



Società Chimica Italiana
Divisione di Spettrometria
di Massa



XXII International Mass Spectrometry Conference

Florence (Italy) - August 26-31, 2018



ABSTRACT BOOK



Società Chimica Italiana
Divisione di Spettrometria
di Massa



XXII International Mass Spectrometry Conference

ABSTRACT BOOK

Florence (Italy) – August 26-31, 2018

TABLE OF CONTENTS

<i>IMSF Executive Committee, IMSF National Affiliates</i>	4
<i>IMSC 2018 Scientific Committee</i>	5
<i>Oral communications</i>	6
<i>a. Sunday, August 26</i>	7
<i>b. Monday, August 27</i>	8
<i>c. Tuesday, August 28</i>	61
<i>d. Wednesday, August 29</i>	114
<i>e. Thursday, August 30</i>	164
<i>f. Friday, August 31</i>	214
<i>Poster communications</i>	240
<i>a. Monday, August 27</i>	241
<i>b. Tuesday, August 28</i>	510
<i>c. Wednesday, August 29</i>	775
<i>d. Thursday, August 30</i>	1019
<i>Author Index</i>	1246

IMSF EXECUTIVE COMMITTEE, IMSF NATIONAL AFFILIATES

International Mass Spectrometry Foundation (IMSF)

Executive Committee

- ✓ Prof. Catherine E. Costello (USA) – President
- ✓ Prof. G. John Langley (UK) – Vice President
- ✓ Prof. Stephan Blanksby (AU) – Secretary
- ✓ Prof. Ron Heeren (NL) – Treasurer
- ✓ Prof. Marcos N. Eberlin (BR) – Past President
- ✓ Prof. Gianluca Giorgi (IT) – Vice President (Conference)
- ✓ Prof. Marcos N. Eberlin (BR) – Vice President (Education)
- ✓ Prof. Michal Sharon (IL) – Education Officer
- ✓ Prof. Gianluca Giorgi (IT) – Region A Representative
- ✓ Prof. Myeong Hee Moon (KS) – Region B Representative
- ✓ Prof. Scott A. McLuckey (USA) – Region C Representative
- ✓ Prof. Marcos N. Eberlin (BR) – Region D Representative

National Affiliates

- | | |
|--------------------------------------|--|
| ✓ Argentina – Rosa Erra Balsells | ✓ Korea – Yong-Hyeon Yim |
| ✓ Australia/New Zealand – Gavin Reid | ✓ Netherlands – Manfred Wuhrer |
| ✓ Austria – Günter Allmaier | ✓ Norway – Einar Uggerud |
| ✓ Belgium – Jeff Rozenski | ✓ Peoples Republic of China – Jin-Ying Li |
| ✓ Brazil – Monica Padilha | ✓ Poland – Piotr Stefanowicz |
| ✓ Canada – Lekha Sleno | ✓ Portugal – Maria Helena Florencio |
| ✓ Croatia – Mario Cindric | ✓ Romania – Zaharie Moldovan |
| ✓ Czech Republic – Jan Preisler | ✓ Russia – Albert Lebedev |
| ✓ Denmark – Steen Pontoppidan | ✓ Serbia & Montenegro – Mila Lausevic |
| ✓ Egypt – Ezzat Selim | ✓ Singapore – Lin Qingsong |
| ✓ Finland – Martin Söderström | ✓ Slovakia – Vladimir Patoprsty |
| ✓ France – Aura Tintaru | ✓ Slovenia – Helena Prosen |
| ✓ Germany – Andrea Sinz | ✓ South Africa – Egmont Rohwer |
| ✓ Hellenic – Despina Tsipi | ✓ Spain – Encarnación Moyano |
| ✓ Hong Kong – April S.Y. Wong | ✓ Sweden – Jonas Bergquist |
| ✓ Hungary – Laszlo Lelik | ✓ Switzerland – Bertran Gerrits |
| ✓ India – Suresh Aggarwal | ✓ Taiwan – Pao-Chi Liao |
| ✓ Ireland – Philip White | ✓ Ukraine – Marina Kosevich |
| ✓ Israel – Dana Reichmann | ✓ United Kingdom – Gavin O'Connor |
| ✓ Italy – Donatella Caruso | ✓ United States of America – Jennifer Brodbelt |
| ✓ Japan – Mitsuo Takayama | |

22nd IMSC SCIENTIFIC COMMITTEE

Chair: Gianluca Giorgi University of Siena

Scientific Committee:

Giuseppe Avellone	University of Palermo	Roberta Galarini	IZS Umbria-Marche, Perugia
Gianluca Bartolucci	University of Florence	Giancarlo la Marca	Meyer Hospital, Florence
Giuliana Bianco	University of Basilicata, Potenza	Fulvio Magni	University Milano Bicocca, Milan
Donatella Caruso	University of Milan	Giorgio G. Mellerio	University of Pavia
Leopoldo Ceraulo	University of Palermo	Gloriano Moneti	University of Florence
Francesco De Angelis	University of L'Aquila	Giovanni Sindona	University of Calabria
Simonetta Fornarini	Sapienza, University of Rome	Michele Suman	Barilla, Parma
		Pietro Traldi	CNR, Padova



Gianluca Giorgi



Giuseppe Avellone



Gianluca Bartolucci



Giuliana Bianco



Donatella Caruso




Leopoldo Ceraulo



Francesco De Angelis



Simonetta Fornarini

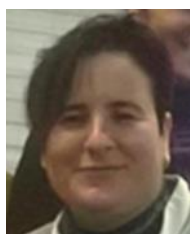

Roberta Galarini



Giancarlo la Marca



Fulvio Magni



Giorgio G. Mellerio



Gloriano Moneti



Giovanni Sindona



Michele Suman



Pietro Traldi

ORAL COMMUNICATIONS

PL01: Plenary Lecture

PL01 **AT THE INTERSECTION BETWEEN CHEMISTRY AND ART: SCIENTIFIC RESEARCH FOR THE STUDY AND PRESERVATION OF CULTURAL HERITAGE****Marco Leona***Department of Scientific Research, The Metropolitan Museum of Art, 1000 Fifth Avenue, New York, NY 10028***Keywords:** cultural heritage; art; Raman; x-ray fluorescence; mass spectrometry

The creation of a work of art is a process often depicted as untouched by material constraints. Pure, unfettered inspiration, without the limitations imposed by physics. It should be evident, however, that only with materials having defined physical properties a certain esthetic effect can be obtained. Works of art are complex technological objects: their appearance and their performance over time are the results of a careful selection of materials and techniques. Thus, artists' choices can be interpreted not just in the classical framework of art history, but in a holistic approach encompassing esthetics and technology. In absence of direct testimony from their often anonymous creators, masterpieces of the past become key material evidence in our efforts to retrace the working practices of artists, the technical knowledge and esthetic preferences of ancient cultures, and the complex trade patterns linking distant empires.

The quest for color, for beauty is a fundamental constant in human history, as well as an economic engine, and an inspiration for art and technology alike. While we easily recognize that the desire to produce richer and more lively images is a clear driver to seek new materials, on the other side of the world, in mines or forests, or in the crucibles of the alchemist (or of his modern counterpart, the chemist), we often overlook how new materials lead to new esthetics.

Scientific investigations have expanded significantly our understanding of works of art, often upending traditional interpretations. A typical example is the rediscovery and investigation of the original polychromy of ancient Greek and Roman sculpture, once a contested topic and now an established research field. Many more cases highlight the importance of the work of cultural heritage scientists, not only in investigating works of art, but also in developing the techniques that make these studies possible. Examples include the discovery of the earliest organic pigment, in a 4,000 years old Egyptian leather fragment; new insights into the global trade of dyes in the Middle Ages; and closer to us, the introduction of synthetic pigments and dyes in the 19th centuries, from the craze for mauveine to Van Gogh's unfortunate passion for geranium lakes, passing by Hokusai's experiments with Prussian blue in the *Thirty-six views of Mount Fuji*, and the color revolution prompted by the arrival into Japan of the aniline dyes in the Meiji period.

Using advanced analytical chemistry, from x-ray fluorescence spectrometry to Raman and surface-enhanced Raman spectroscopy, from electron microscopy to mass spectrometry, we can build a *material history of art* based on physical evidence gathered from masterpieces spanning four millennia, from ancient Egypt to Modernity, demonstrating that technological developments were not only readily embraced, and often prompted by artists and their audiences, but also that they in turn created new forms of expression, in a continuous dialogue between materials and meanings.

PL02: Plenary Lecture

Chair: Ron M. A. Heeren

PL02 ADVANCES IN MALDI IMAGING MASS SPECTROMETRY: MOLECULAR MICROSCOPY IN THE NEW AGE OF BIOLOGY AND MEDICINE

Richard Caprioli

Depts. of Biochemistry, Chemistry, Pharmacology and Medicine Vanderbilt University School of Medicine

MALDI Imaging Mass Spectrometry (IMS) has been employed in studies of a wide variety of biologically and medically relevant research projects. This technology has seen significant advances in both the instrumentation and in the fields of basic biological investigations and medically relevant applications including patient biopsies. Several of these applications will be presented in this presentation. In addition, new technological advances will be described, including sample preparation and instrumental performance to achieve images at high spatial resolution (1-5 microns) and at high speeds in order to minimize image acquisition time. Instrumentation improvements include those developed for both MALDI FTICR MS and MALDI TOF mass spectrometers and will be demonstrated by specific applications. Finally, new biocomputational approaches will be discussed that deals with the high data dimensionality of IMS, image registration with microscopy, and our implementation of 'image fusion' in terms of predictive integration of MS images with other imaging modalities.

MOr.01 - Nucleic acids

Chairs: Daniele Fabris, Giuliana Bianco

MOr.01-1 Keynote: STRUCTURES, ENERGETICS AND REACTIVITY OF NUCLEIC ACIDS AND THEIR MOLECULAR BUILDING BLOCKS VIA TANDEM MS AND COMPUTATIONS APPROACHES

Mary T. Rodgers

Wayne State University, Detroit, USA

Keywords: Collision-induced dissociation; Glycosidic bond stability; Infrared multiple photon dissociation action spectroscopy; Nucleic acids; Nucleosides

Introduction

Synergistic tandem mass spectrometry and computational chemistry approaches enable questions related to the structure, energetics and reactivity of nucleic acids and their building blocks to be investigated. Here the effects of modification and the local environment on nucleoside conformation and glycosidic bond stability and the reactivity of nucleic acids with potential drug candidates (amino glycoside, peptide and metal-ligand complexes) are presented.

Methods

Several types of tandem mass spectrometers (GIBMS, QITMS, and FT-ICR MS) and approaches including: threshold, energy-resolved, and sustained off-resonance collision-induced dissociation (TCID, ER-CID, and SORI-CID) and wavelength-resolved infrared multiple photon dissociation

(IRMPD) action spectroscopy, complemented by electronic structure calculations of the stable structures, energetics, IR spectra, and reactions mechanisms are employed.

Results

Combined IRMPD action spectroscopy, TCID, ER-CID, and computational chemistry studies of the protonated and sodium cationized forms of the canonical DNA and RNA nucleosides as well as a variety of naturally occurring and synthetically modified nucleosides have been performed. The preferred sites of protonation or sodium cationization, nucleobase orientations, sugar puckerings, conformers populated by electrospray ionization, and glycosidic bond stabilities are determined. The reactivity and preferred sites of binding of potential drug candidates with pre-micro and ribosomal RNA hairpins are elucidated using top-down sequencing and SORI-CID approaches. The RNAs investigated include: Helix 69, a hairpin RNA structure containing 19 nucleotides in the large subunit of the ribosome, and TAR, the HIV trans-activation response element, a hairpin RNA that acts as the binding site for the Tat protein. Drug candidates examined include: cisplatin and a series of amino acid-linked cisplatin derivatives, amino glycosides, and peptides.

Conclusions

The structures and glycosidic bond stabilities of protonated and sodium cationized canonical and modified nucleosides provide insight into structure-function of these genetics building blocks, evolutionary aspects of roles of modification, and impact of unnatural modifications. The binding and reactivity of drugs to RNAs provide insight into structure-function of functional RNAs, and sample preparation lessons for achieving quality mass spectral analyses.

Novel Aspect

Effects of local environment/modifications on structures/stabilities of nucleosides are elucidated. Reactivity/binding sites of drug candidates to RNAs are determined.

MOr.01-2 **TANDEM MASS SPECTROMETRIC SEQUENCING OF ISOMERIC METHYLNUCLEOSIDES ON RNA TOWARD THE UNDERSTANDING OF EPITRANSCRIPTOME**

Hiroshi Nakayama⁽¹⁾ - ***Yoshio Yamauchi***⁽²⁾ - ***Toshiaki Isobe***⁽³⁾ - ***Masato Taoka***⁽³⁾

⁽¹⁾ *Biomolecular Characterization Unit, RIKEN Center for Sustainable Resource Science, Wako* - ⁽²⁾ *Department of Chemistry, Tokyo Metropolitan University, Hachioji* - ⁽³⁾ *Department of Chemistry, Tokyo Metropolitan University, Hachioji*

Keywords: RNA, post-transcriptional modification, epitranscriptome, tandem mass spectrometry, methylation

Introduction

The discovery of dynamic post-transcriptional modifications (PTMs) of ribonucleic acid (RNA) including N6-methyladenosines affects transcriptome, emerging field of 'epitranscriptomics'. Although its analysis requires a general sequencing method of modified nucleosides on RNA, conventional methods fail to provide such information. We present here an MS3-based method to sequence positional isomers of methyl nucleosides in RNAs.

Methods

Base-methyl adenosines (mAs), cytidines (mCs), guanosines (mGs), and uridines (mUs) were analyzed by pseudo MS3 [1] with a Q-Exactive; a base anion generated by in-source collision-induced dissociation (sCID) of each nucleoside was fragmented by MS/MS to obtain its characteristic pattern. An RNase T1 digest of tRNA-Phe was sequenced by nano-flow liquid

chromatography-MS/MS in data-dependent [2] and pseudo MS3 [1] modes. The data were processed by Ariadne [3].

Results

To sequence monomethyl nucleoside isomers in RNA, we extended the pseudo MS3 method [1]. We inspected sCID spectra of RNase-digested oligonucleotides containing mono-methyl nucleoside in detail and confirmed an ion at m/z 225.02 exclusively from 2'-O-methylated oligonucleotides and another ion at m/z 148.07 specific for the methyladenine (mAde)-containing oligonucleotides. To further distinguish those mAde isomers, we examined their fragmentations using commercially available mAs; we found each CID generated mAde anion exhibited a characteristic product-ion pattern for each mA isomer. Likewise, we defined the MS3-derived signature patterns for mCs, mGs and mUs. We then validated the method by determining the chemical structure of yeast tRNA-Phe, which is mono-methylated at 8 positions out of 76 nucleotides. The MS3-based analysis of its RNase T1 digest (500 fmol) unambiguously confirmed the site and positional isomerism of all monomethylnucleosides in the tRNA. The method will be applied to sequence methylnucleosides on any RNAs.

Conclusions

We have established signature fragmentation patterns for mAs, mCs, mGs and mUs and, using them in the pseudo MS3 method, sequenced monomethylnucleosides on RNA. In addition, the method should extend to the other nucleosides, e. g. dimethyl ones. Because the method is applicable to a sub-picomole quantity of RNA sample, it will serve as a useful tool for analyzing the structure/function of cellular RNAs and toward the comprehensive understanding of epitranscriptome.

Novel Aspect

A pseudo MS3 analysis allows discrimination of a monomethylnucleoside from its positional isomers with a sub-picomole quantity of RNA.

References

1. Yamauchi Y., Nobe Y., Izumikawa K., Higo D., Yamagishi Y., Takahashi N., Nakayama H., Isobe T., Taoka M., *Nucleic Acids Res*, 44, e59(2016)
2. Nakayama H., Yamauchi Y., Taoka M., Isobe T., *Anal Chem*, 87, 2884-91(2015)
3. Nakayama H., Akiyama M., Taoka M., Yamauchi Y., Nobe Y., Ishikawa H., Takahashi N., Isobe T., *Nucleic Acids Res*, 37, e47 (2009).

MOr.01-3 TEAD1: YIN AND YANG OF THE CELL REGULATION

Petr Novak⁽¹⁾ - **Ruzena Liskova**⁽¹⁾ - **Lukas Slavata**⁽¹⁾ - **Jiri Cerny**⁽²⁾ - **Karel Valis**⁽³⁾

⁽¹⁾ Institute of Microbiology, Laboratory of Structural Biology, Prague - ⁽²⁾ Institute of Biotechnology, Laboratory of Biomolecular Recognition, Vestec - ⁽³⁾ Institute of Microbiology, Laboratory of Structural Biology, Vestec

Keywords: structural proteomics, protein/DNA interaction, chemical cross-linking, H/D exchange, chromatin immunoprecipitation.

Introduction

TEAD transcription factors mediate gene expression regulation through interactions with their DNA response M-CAT motif. They are active during growth and development and activate gene expression of proteins responsible for cell proliferation, differentiation or apoptosis. To understand the principles of interaction between TEAD1 and M-CAT motifs originating from regulatory regions of human genes, several structural mass spectrometry techniques were used.

Methods

TEAD1 DNA binding domain (TEAD1-DBD) was expressed in *E.coli*. Fluorescence anisotropy was used to determine K_D for TEAD1-DBD/M-CAT complexes. H/D exchange was carried out for DNA-free TEAD1-DBD or in complex with each M-CAT motif. DSA in non- ^{13}C -labelled forms was used for quantitative cross-linking. LCMS analysis was performed on 15T solarix. Structural models were built using MODELLER and *in vivo* quantification of M-CAT occupancy was determined byChIP-Q-PCR.

Results

According to K_D values, tested M-CAT motifs can be divided into two groups one with ten times higher affinity to TEAD1-DBD than the other. H/D exchange results revealed differences in deuterium uptake in helix H3, part of helix H2 and in the loop connecting them, identifying this as DNA-binding region. In correlation with dissociation constants, protection from deuteration in higher-affinity complexes was more intense than in low-affinity complexes. Quantitative chemical cross-linking resulted in 14 distance restrains. Three of them were not affected by presence of DNA, two were favored in presence of M-CAT and eight constrains were discriminated by DNA binding. Cross-link formation ability of lysines K57 and K88 significantly decreased in complexed state. Again, the effect of high-affinity duplexes was more significant than the effect of low-affinity ligands. H/D exchange and cross-linking results allowed creation of TEAD-DBD/M-CAT structural models. The chromatin immunoprecipitation confirmed such observation in living cells.

Conclusions

Experimental data uncovers higher affinity of TEAD1 to 5'-3' than 3'-5' oriented M-CAT motifs and enables molecular docking of TEAD1/M-CAT complex. Based on structural models TEAD1 recognizes its response motif by the shape of DNA major groove and the complex is stabilized by specific aminoacid-nucleotide interactions. TEAD1 binds to inverted M-CATs 180° degrees rotated, less and weaker bonds are formed. Such observation correlates with the *in vivo* situation.

Novel Aspect

Structural mass spectrometry sheds a light on the tumorigenesis orchestrated by TEAD1 transcriptional factor.

For information please contact: pnovak@biomed.cas.cz

MOr.01-4 **HSA-MIR-1587 G-QUADRUPLEX FORMATION AND DIMERIZATION INDUCED BY NH_4^+ , MOLECULAR CROWDING ENVIRONMENT AND JATRORRHIZINE DERIVATIVES**

Wei Tan⁽¹⁾ - **Long Yi**⁽²⁾ - **Zhentao Zhu**⁽²⁾ - **Lulu Zhang**⁽¹⁾ - **Jiang Zhou**⁽¹⁾ - **Gu Yuan**⁽¹⁾

⁽¹⁾ Peking University, Chemistry, Beijing - ⁽²⁾ Beijing University of Chemical Engineering, Beijing Key Laboratory of Bioprocess, Beijing

Keywords: miR-1587; G-quadruplex; dimerization; jatrorrhizine derivatives; ESI-MS

Introduction

G-quadruplexes widely exist in genome and transcriptome [1], which have been considered as promising therapeutic targets [2]. However, reports about RNA G-quadruplexes are relatively fewer [3,4] and factors that influence high-order RNA G-quadruplex formation are still under to be investigated. In this study, ESI-MS was utilized to study the formation and properties of a RNA G-quadruplex formed by miR-1587.

Methods

Most of the ESI-MS experiments were performed on a Finnigan LCQ Deca XP Plus ion-trap mass spectrometer in the negative ion mode. The high resolution ESI-MS spectra were collected on a Bruker Solarix XR FT-ICR mass spectrometer with a 9.4T magnet. All the CD experiments were carried out on a JASCOJ-815 CD spectrometer equipped with a Peltier junction temperature controlled cell holder and a 1.0 cm cell.

Results

A guanine-rich human mature microRNA, miR-1587, was discovered to form stable intramolecular G-quadruplexes in the presence of K^+ , Na^+ and low concentration of NH_4^+ (25 mM) by electrospray ionization mass spectrometry (ESI-MS) combined with circular dichroism (CD) spectroscopy. Furthermore, under high concentration of NH_4^+ (100 mM) or molecular crowding environment, miR-1587 formed a dimeric G-quadruplex through 3'-to-3' stacking of two monomeric G-quadruplex subunits with one ammonium ion sandwiched between the interfaces. In contrast, DNA-1587 and dU-DNA-1587 sequences also formed three-layer G-quadruplexes, but was less stable and could not form dimeric G-quadruplex structures. Specifically, two synthesized jatrorrhizine derivatives with terminal amine groups could also induce the dimerization of miR-1587 G-quadruplex and formed 1:1 and 2:1 complexes with the dimeric G-quadruplex. In contrast, jatrorrhizine could bind with the dimeric miR-1587 G-quadruplex, but could not induce dimerization of miR-1587 G-quadruplex.

Conclusions

In this study, the formation and property of intramolecular miR-1587 G-quadruplex were investigated. High concentration of NH_4^+ , molecular crowding environment and two jatrorrhizine derivatives were discovered to be capable of inducing the formation of dimeric miR-1587 G-quadruplex. The results expanded our knowledge of RNA G-quadruplexes and provided new insights into the targeting of higher-order structures of miRNAs to regulate their functions.

Novel Aspect

A dimeric miRNA G-quadruplex was first found by ESI-MS under high concentration of NH_4^+ or molecular crowding environments.

References

1. Burge S., Parkinson G. N., Hazel P., Todd A. K., Neidle S., *Nucleic Acids Res*, 34, 5402 (2006).
2. Collie G. W., Parkinson G. N., *Chem Soc Rev*, 40, 5867 (2011).
3. Collie G. W., Haider S. M., Neidle S., Parkinson G. N., *Nucleic Acids Res*, 38, 5569 (2010).
4. Martadinata H., Phan A. T., *Biochemistry*, 53, 1595 (2014).

MOr.01-5 BIOPHYSICAL CHARACTERIZATION OF APTAMER-METABOLITE INTERACTIONS USING ION MOBILITY MASS SPECTROMETRY

Stefano Piccolo⁽¹⁾ - **Elise Daems**⁽²⁾ - **Frank Sobott**^(2,3,4) - **Valérie Gabelica**⁽¹⁾

⁽¹⁾ Laboratoire ARNA, INSERM U1212, CNRS UMR 5320, Univ. Bordeaux, IECB, ARNA Laboratory, F-33600 Pessac, France - ⁽²⁾ Biomolecular and Analytical Mass Spectrometry, University of Antwerp, 2020 Antwerp, Belgium - ⁽³⁾ Astbury Centre for Structural Molecular Biology, University of Leeds, Leeds, United Kingdom - ⁽⁴⁾ School of Molecular and Cellular Biology, University of Leeds, Leeds, United Kingdom

Keywords: Native MS, Aptamers, Nucleic Acids, CIU, Ion Mobility

Introduction

Aptamers are Nucleic Acids (NAs) capable to bind selectively to a ligand, or to a family of molecules. In aptamer-based technology (e.g., drug delivery or sensors), it is crucial to understand how binding is occurring, to quantify affinities, and how they can be tuned. We

interrogated several RNA and DNA aptamers by native MS and ion mobility, to understand how they fold upon binding and which quantifications limits are attainable.

Methods

All NAs are analysed in MS-compatible solutions such as TMAA or NH₄OAc with Mg²⁺ and K⁺(when needed) to mimic bio-relevant conditions. An Agilent 6560 ESI-IMS-Q-TOF in negative mode is used for most experiments and CCS were determined in helium. A Waters Synapt G2 HDMS is used for IMS comparisons on cocaine aptamers. A MicroCal ITC 200 is used as an independent measurement of K_Ds. UV-melting and CD in solution is used to get information about NAs folding state in solution.

Results

Among the different systems tested, only DNA Cocaine aptamers MN4, MN19 and OR8 [1] and RNA Tetracycline aptamer [2] are giving significant conformational change upon binding by IM-MS. Although for tetracycline complex and free-aptamer two distinct conformations are present, cocaine aptamers are showing a transition between two conformational ensembles only for the free-aptamer ion. For cocaine aptamers CIU on charge state 7- is best revealing their relative stabilities in gas phase in function of length.

Unfortunately, UV-Tm and CD do not allows to predict the conformational space in gas phase for the NAs tested. But we had recently evidence that collapsing can occur on duplex structures in gas phase [3].

Besides their shape, aptamer-ligand affinities are determined by MS-titration using a dT₆ internal standard [4], and we defined the lower and higher K_D attainable. ITC and MS K_D values are comparable, both in TMAA and NH₄OAc. For Malachite green and Cocaine aptamer MS data were precious to clarify the number of sites on ITC.

Conclusions

Hairpin patterns in Aptamers are prone to compact and lose the in solution structure “memory” during ESI, if the charge states are too low. Here we show for the first time two cases, Cocaine (bound-form) and Tetracycline aptamers, where bound and unbound-form are giving distinct IMS signature at some charge states. We also suggest native MS as an independent quantification technique and to reconcile puzzling data from ITC and NMR.

Novel Aspect

We discuss the advantages and limitations of native MS and ion mobility of DNA and RNA aptamers and their ligand complexes

References

1. Neves et al., *ACS Sens.* **2017**, 2, 1539
2. Reuss et al., *Biophys. J.* **2014**, 107, 2962.
3. Porrini et al., *ACS Cent. Sci.* **2017**, 3, 454.
4. Gabelica et al., *Anal. Chem.* **2009**, 81, 6708.

This work is funded by the EU H2020 (Marie Curie ITN “MetaRNA”)

MOr.02 - Clinical proteomics

Chairs: Giuseppina Maccarrone, Donatella Caruso

MOr.02-1 Keynote: IMMUNOASSAYS IN MULTIPLEX FOR BIOMARKER DISCOVERY AND VALIDATION

Thomas O. Joos

University of Tübingen Reutlingen

Keywords: Biomarker, Immunoassays, Safety biomarkers, Signaling

Introduction

Array-based assay systems allow the analysis of hundreds of molecular parameters in a single experiment. Within the last decade multiplexed protein analytical assays achieved robust analytical performance and enable to screen for a multitude of parameters using minimal amounts of sample material. There are a variety of technologies and methods available to address assay requirements in terms of covering analyte concentration, sample variety and sample origin.

Methods

Two assay formats addressing unmet needs in biomarker discovery developed at NMI will be presented.

Our DigiWest assay format transfers the western blot to a bead-based microarray platform extending traditional western blotting to enable the rapid quantification of hundreds of specific proteins and protein modifications in complex samples (1).

The second assay format combines sensitive and selective immunoprecipitation coupled with mass spectrometric readout. The use of peptide group-specific antibodies called MS-based immunoassays employing such terminal specific antibodies are established for several applications, including a number of cross-species drug induced organ injury assays that represent truly translational biomarker assays.(2)

Results

The DigiWest high-throughput western blot approach allowed us to identify alterations in cellular signal transduction that occur during the development of resistance to the kinase inhibitor Lapatinib. A total of 185 western blot equivalents were performed. Differential phosphorylation was observed for a variety of proteins; a substantial increase in phosphorylated Aurora A, p27 and IRAK4 was detected; for other proteins a decrease in phosphorylation was observed (e.g. MEK1, Src, Erk1/2, different sites in EGFR). (1)

For our immunoaffinity mass spectrometry assays the quantification of drug induced kidney biomarker candidates in canine and monkey urine samples will be presented. The expression profiles of osteopontin, retinol binding protein 4, clusterin and KIM-1 in urine samples correlates with histopathology and conventional clinical pathology parameters when compared to tissue samples from the same animal. The biomarker assays are more sensitive than blood urea nitrogen (BUN) and serum creatinine (sCr) to detect early drug induced kidney injury. (4)

Conclusions

Biomarker discovery and subsequent validation require robust technologies and assays combined with an excellent study design to include appropriate sample sets. This is a prerequisite for any biomarker qualification process. Our new assay formats enable the generation of such relevant data sets and can be implemented into a generic qualification strategy. Such new biomarkers will allow early identification, assessment and management of drug-induced injuries throughout R&D as well as biomarkers relevant for precision medicine.

References

1. Treindl F, Ruprecht B, Beiter Y, Schultz S, Döttinger A, Staebler A, Joos TO, Kling S, Poetz O, Fehm T, Neubauer H, Kuster B, and Templin MF. *Nat Commun.* 2016; 7: 12852.
2. Weiß F, Schnabel A, Planatscher H, van den Berg BH, Serschnitzki B, Nuessler AK, Thasler WE, Weiss TS, Reuss M, Stoll D, Templin MF, Joos TO, Marcus K, Poetz O. *Sci Rep.* 2015; 5:8759
3. Weiss F, van den Berg BH, Planatscher H, Pynn CJ, Joos TO, Poetz O. *Biochimica et Biophysica Acta* 2013.
4. Gautier JC, Zhou X, Yang Y, Gury T, Qu Z, Palazzi X, Léonard JF, Slaoui M, Veeranagouda Y,

Guizon I, Boitier E, Filali-Ansary A, van den Berg BH, Poetz O, Joos T, Zhang T, Wang J, Detilleux P, Li B. *ToxicolApplPharmacol.* 2016 Jul 15;303:1-10

MOr.02-2 **SHEDDING NEW LIGHT ON SPINAL CORD INJURY**

Michel Salzet⁽¹⁾ - **Jusal Quanico**⁽²⁾ - **Melodie Anne Karnoub**⁽²⁾ - **Celine Meriaux**⁽²⁾ - **Zahra Laouby**⁽²⁾ - **Franck Rodet**⁽³⁾ - **Dasa Cizkova**⁽⁴⁾ - **Isabelle Fournier**⁽²⁾ - **Richard Assaker**⁽²⁾

⁽¹⁾ U1192 Inserm, PRISM Laboratory, Univeristy Lille, Villeneuve d'Ascq - ⁽²⁾ U1192 Inserm, PRISM Laboratory, University Lille, Villeneuve d'Ascq - ⁽³⁾ U1192 Inserm, PRISM Laboratory, University Lille, Villeneuve D'Ascq - ⁽⁴⁾ Institute of Neuroimmunology, Slovak Academy of Sciences, Bratislava

Keywords: Spinal Cord Injury, Spatio temporal Proteomic, Tempotally 3D MALDI MSI, Antibody Dependant Neuritogenesis Modulation, Stem Cells

Introduction

Spinal cord injury (SCI) belongs to currently incurable disorders of the CNS and is accompanied by permanent health consequences-disability. In order to mimic a SCI, a balloon-compressive technique was used at thoracic Th8-9 spinal level in adult rat.

Methods

4D MALDI Imaging, lipidomics, Saptiotempral tissue microproteomics were undertaken combined with confocal imaging. Exosomes from Stem Cells, functionalized biomaterial, Rho A inhibitor have been tested in pre-clinical way to develop a smart biomaterial.

Results

We determined the spatio-temporal events occurring in acute phase after SCI. Caudal segment has clearly been detected as the therapeutic target [1,2]. We then assessed in a rat SCI model the in vivo impact of a sustained RhoA inhibitor administered in situ via functionalized-alginate scaffold [3,4] . In order to decipher the underlying molecular mechanisms involved in such a process, an in vitro neuroproteomic-systems biology platform was developed in which the pan-proteomic profile of the dorsal root ganglia (DRG) cell line ND7/23 DRG was assessed in a large array of culture conditions using RhoAi and/or conditioned media obtained from SCI ex-vivo derived spinal cord slices. A fine mapping of the spatio-temporal molecular events of the RhoAi treatment in SCI was performed. The data obtained allow a better understanding of regeneration induced above and below the lesion site.

Conclusions

Results notably showed a time-dependent alteration of the transcription factors profile along with the synthesis of growth cone-related factors (receptors, ligands, and signaling pathways) in RhoAi treated DRG cells and involvement of IgG by binding to their receptors on the DRG cells

Novel Aspect

We established a novel origin of IgG, their role in neurites outgrowth modulation, and developed a smart biomaterial for treating SCI.

References

1. Cizkova, D., Le Marrec-Croq, F., Franck, J., Slovinska, L., Grulova, I., Devaux, S., Lefebvre, C. Fournier, I. Salzet, M. *Frontiers in Cell. Neurosc.* 8:105 (2014)
2. Devaux S., Cizkova D., Qunaico J., Franck J., Ntaf S., Pays L., Hauberg-Lotte L., Mass P., Hobart J.H., Kobeissy F., Merieux, C., Wisztorski M., Slovinska L., Blaska J., Cifankova V., Fournier I., Salzet M. *Mol Cell Proteomics.* (15(8):2641-70 (2016)
3. Devaux S, Cizkova D, Mallah K, Karnoub MA, Laouby Z, Kobeissy F, Blasko J, Nataf S, Pays L,

- Mériaux C, Fournier I, Salzet M. *Mol Cell Proteomics*. 16(8):1394-1415 (2017)
4. Cizkova, D., Devaux, S., Le Marrec-Croq, F., Franck, J., Slovinska, L., Rosocha, J, Spakova, T., Lefebvre, C. Fournier, I. Salzet, M. Modulatory properties of mesenchymal stem cells on BV2 microglia after spinal cord injury secretome activation: Proteomic and physiological studies. *Sci. Report*. (2014) 4:7514.
 5. Grulova I, Slovinska L, Blaško J, Devaux S, Wisztorski M, Salzet M, Fournier I, Kryukov O, Cohen S, Cizkova D. *Sci Rep*. 5:13702 (2015)
 6. Cizkova D, Cubinkova V, Smolek T, Murgoci AN, Danko J, Vdoviakova K, Humenik F, Cizek M, Quanico J, Fournier I, Salzet M. *Int J Mol Sci*. 9(3). pii: E870 (2018)

MOr.02-4 IDENTIFICATION OF PROTEOMIC AND LIPIDOMIC MARKERS IN INFLAMMATORY BOWEL DISEASES BY MS-BASED TECHNIQUES AND MULTIVARIATE STATISTICAL ANALYSIS

Elisa Robotti⁽¹⁾ - ***Marcello Manfredi***⁽¹⁾ - ***Eleonora Conte***⁽²⁾ - ***Elettra Barberis***⁽¹⁾ - ***Arianna Buzzi***⁽¹⁾ - ***Valeria Caneparo***⁽³⁾ - ***Daniela Cecconi***⁽⁴⁾ - ***Jessica Brandi***⁽⁴⁾ - ***Ester Vanni***⁽⁵⁾ - ***Santo Landolfo***⁽⁶⁾ - ***Marisa Gariglio***⁽⁷⁾ - ***Emilio Marengo***⁽¹⁾ - ***Marco De Andrea***⁽⁷⁾

⁽¹⁾ *University of Piemonte Orientale, Department of Sciences and Technological Innovation, Alessandria* - ⁽²⁾ *Isalit Srl, Novara* - ⁽³⁾ *Medical School of Novara, Novara* - ⁽⁴⁾ *University of Verona, Department of Biotechnology, Verona* - ⁽⁵⁾ *University of Torino, Department of Medical Sciences, Torino* - ⁽⁶⁾ *Medical School of Turin, Department of Public Health and Microbiology, Torino* - ⁽⁷⁾ *University of Piemonte Orientale, Dimet, Novara*

Keywords: Inflammatory bowel diseases; proteomics; lipidomics; multivariate statistical analysis; biomarkers identification

Introduction

The causes of Inflammatory Bowel Diseases (IBDs) are still little known, the prevailing hypothesis being an immunologic reaction towards antigens. The diagnosis of IBDs is usually carried out by quite invasive instrumental techniques; some plasma biomarkers have been proposed but much can be done to clarify the etiology of the pathology and identify reliable biomarkers with predictive ability.

Methods

Plasma proteins and fatty acids from controls and patients affected by Chron's disease (CD) and Ulcerative Colitis (UC) were analysed using mass spectrometry-based techniques. Low-abundant plasma proteins were quantified after depletion of the 14 most abundant and tryptic digestion, followed by uHPLC-QTOF-MS/MS. Several fatty acids were simultaneously analysed by Gas Chromatography Mass Spectrometry (GC-MS), after derivatization.

Results

uHPLC-QTOF-MS/MS was applied exploiting the SWATH-MS acquisition, a high throughput label-free method for protein quantitation combining shotgun proteomics with the quantitative accuracy and reproducibility of selected reaction monitoring (SRM). The data collected were further processed by multivariate statistical tools to identify markers of CD and UC but also identify markers able to discriminate UC vs CD samples. Partial least squares-discriminant analysis (PLS-DA) coupled to a variable selection procedure in cross-validation was applied to this purpose, in order to identify the most relevant panel of markers with the best predictive ability. Proteomic and lipidomic data were also integrated to provide new combined biomarkers including both proteomic and lipidomic markers. The analysis allowed the identification of panels of markers with a very good predictive ability and pointed out the importance of including fatty acids.

Conclusions

The combination of proteomic and lipidomic data allowed to obtain panels of biomarkers with

improved predictive ability, pointing out the importance of fatty acids, combined with proteins, in the identification of IBDs. Bioinformatic tools applied to the panels of identified markers allowed the identification of deranged pathways (up and/or down regulated).

Novel Aspect

The combination of proteomic and lipidomic data provided panels of markers with improved predictive ability with respect to the only proteomic analysis.

MOr.02-5 A NOVEL MASS SPECTROMETRY BASED ASSAY FOR THE DIAGNOSIS AND TYPING OF SYSTEMIC AMYLOIDOSIS

Paul Boersema⁽¹⁾ - ***Ting Huang***⁽²⁾ - ***Johan Bijzet***⁽³⁾ - ***Christian Zimmerli***⁽⁴⁾ - ***Olga Vitek***⁽²⁾ - ***Bouke Hazenberg***⁽³⁾ - ***Paola Picotti***⁽⁵⁾

⁽¹⁾ *ETH Zurich, Institute of Molecular Systems Biology, Zurich* - ⁽²⁾ *Northeastern University, College of Computer and Information Science, Boston* - ⁽³⁾ *University Medical Center Groningen, Rheumatology, Groningen* - ⁽⁴⁾ *ETH Zurich, Institute of Molecular Systems Biology, Zurich* - ⁽⁵⁾ *ETH Zurich, Institute of Molecular and Systems Biology, Zurich*

Keywords: systemic amyloidosis, SRM, diagnosis, fat tissue aspirates, clinical

Introduction

In systemic amyloidoses (SA) different precursor proteins misfold, aggregate as amyloid fibrils and are deposited in various tissues. Correct identification of the precursor amyloid protein is key for successful clinical treatment. However, current diagnostic approaches are lacking. We developed an LC-MS method that assays in one analysis the major specific systemic amyloidosis proteins, their mutations and additional proteins as generic indicators of SA.

Methods

We developed an SRM-MS assay targeting 13 major SA proteins and often-occurring mutations plus 20 additional proteins. To improve peptide detectability and allow for intra-sample and batch normalization, we spiked in heavy isotope labeled synthetic peptide analogs for all targeted peptides and iRT peptides. The final assay, containing 289 precursors and 1134 transitions, was condensed to a single, 35 minute LC-SRM-MS analysis method.

Results

A 1- μ g fraction of the trypsin digested tissue is enough to perform the developed SRM-assay. After target peptide and transition selection and optimization together with retention time scheduling, we could quantitatively measure more than 1000 transitions in a 35 minute LC gradient.

Three sets of in total 226 samples consisting of control samples and patient samples with varying levels of the SA types ATTR, AA, AL-L and AL-K were measured with the SRM-assay. One set was used to create a prediction model based on amyloid protein intensities. The two other sets were used to test and validate the prediction model. The amyloidogenic precursor protein was correctly detected with >80% specificity and sensitivity. The correctness of diagnosis increased with increasing disease severity. Potentially present mutations could be detected simultaneously with the assay.

Conclusions

An amyloid signature of 20 proteins was developed to diagnose the presence of amyloid fibrils in fat tissue aspirates from systemic amyloidosis patients regardless of type and to gauge SA severity.

The high specificity and sensitivity of the assay combined with its relative straightforwardness make it an interesting alternative to the contemporary SA diagnosis pipeline.

Novel Aspect

Single, 35-minute LC-SRM assay simultaneously diagnoses, types and grades 13 different SA types including their potential mutations.

MOr.03 - Ion spectroscopy & Photodissociation

Chairs: Nicolas C. Polfer, Simonetta Fornarini

MOr.03-1 Keynote: NEW FRONTIERS FOR CRYOGENIC CHEMISTRY AND OPTICAL SPECTROSCOPY OF MASS-SELECTED IONS**Mark Johnson**

Department of Chemistry, Yale University, New Haven, CT, 06520, USA

Keywords: Cryogenic ion spectroscopy, reaction intermediates, catalysis

Introduction

Over the past decade, the structural characterization of mass-selected ions with vibrational and electronic spectroscopies has evolved from a niche technique in the hands of a few specialists to a powerful new capability for routine chemical analysis. We present recent developments in this approach that are designed to capture fragile reaction intermediates, determine their structures, and photoinitiate product formation.

Methods

Vibrational and electronic structure of cold ions generated by ambient ionization sources are obtained using cryogenic RF ion traps (octopole, Paul) coupled to a triple focusing time-of-flight photofragmentation mass spectrometer. This hybrid approach is presently being integrated with a Thermo LTQ Orbitrap XL instrument to provide a general secondary analysis capability. Vibrational and electronic spectra of the cold ions are obtained by photoevaporation of weakly bound mass "tags" (typically He, H₂ and N₂) or, in the case of bare ions, by two-color, vibrationally mediated photodissociation.

Results

Two major advantages of the tagging approach are i) high resolution optical spectra can be obtained with the modest power available with table-top laser sources, and ii) the linear spectra of the cold ions can be directly compared with band patterns for candidate structures calculated with widely available electronic structure packages. This scheme has clarified the chemical nature of many molecular ions as well as the docking arrangements of non-covalently bound complexes ranging from polypeptides to unstable reaction intermediates in catalysis. This talk will highlight recent progress in understanding the origin of the strong anharmonicities that are apparent in spectra even when the ions are cooled below 10 K. These effects contribute to broadening and extra features that complicate structural analyses based on fundamentals calculated at the harmonic level. Examples where this occurs typically involve strong hydrogen bonds near charge centers and multi-state Fermi resonances, both of which are common in the vibrational spectra of gas phase peptides. Theoretical markers will be discussed that signal when such complexities are likely to occur. On-going experimental efforts will also be presented involving the development of chemically active ESI ion sources to provide a rational synthetic strategy with which to trap fragile intermediates in the catalytic activation of small molecules. Next-generation opportunities for this rapidly expanding technology will also be noted, such as the prospects for monitoring reaction kinetics in the size-selected solvation regime.

Conclusions

Optical spectra of cold ions yield structures through direct comparison with predictions from electronic structure calculations.

Novel Aspect

Optical spectroscopies of mass selected ions are becoming mainstream capabilities for general chemical analysis.

MOr.03-2 **STRUCTURAL CHARACTERIZATION OF ETD PRODUCT IONS USING IR ION SPECTROSCOPY**

Jos Oomens - Lianne Kempke - Jonathan Martens - Giel Berden

Radboud University, FELIX Laboratory, Nijmegen

Keywords: Electron Transfer Dissociation, ion spectroscopy, peptides, radicals, ion chemistry

Introduction

Peptide sequencing by Electron Transfer Dissociation (ETD) has seen rapid growth especially in top-down applications and for the identification of post-translational modifications (PTM). Understanding its underlying reaction chemistry remains elusive to some extent, in part due the radical-driven dissociation mechanisms and the open-shell product ion structures. Here we apply IR ion spectroscopy to reveal molecular structures of several ETD product ions.

Methods

Experiments are carried out in a modified ETD enabled quadrupole ion trap MS connected to the IR beamline of the free electron laser FELIX. ETD fragment ions of doubly protonated precursor peptides are formed by ion-ion reactions with the fluoranthene radical anion. Individual ETD fragments are mass isolated and spectroscopically investigated with IRMPD spectroscopy [1]. IR spectra are compared with DFT computed spectra to derive structural information.

Results

ETD is known to lead mainly to cleavage of the backbone N-C α bonds in a peptide and results in closed-shell c-type fragments and radical z-type fragments. We have recorded IR spectra for a comprehensive set of c-type ions, mainly to distinguish between possible amide and enol-imine moieties at their C-terminus. Radical z-type ions are also investigated. One of the main questions was whether the radical remains located at the site of cleavage or whether it migrates to a thermodynamically more favorable site. For peptides where ETD occurs at a Cys residue, we observe secondary neutral loss of a thiol radical and spectroscopically establish the formation of a C=C double bond in the so-called w-type fragment ion.

Conclusions

The IR spectra suggest that the c-type ions have an amide moiety at their C-terminal. This may suggest that they are formed via the Utah-Washington mechanism, but post-ETD rearrangement from enol-imine to amide cannot be excluded based on the spectroscopic data. From the spectra of the radical z-ions, we tentatively conclude that radical migration to aromatic side chains does not occur, although other radical migrations have not been exhaustively verified thus far.

Novel Aspect

Systematic investigation of ETD fragment ion structures by IRMPD spectroscopy.

References

1. Martens, J., Grzetic, J., Berden, G., Oomens J., *Nature Commun.* 7, 11754 (2016).

MOr.03-3 STRUCTURES AND PHOTODISSOCIATION ACTION SPECTROSCOPY OF NUCLEIC ACID CATION RADICALS

Frantisek Turecek

University of Washington, Chemistry, Seattle

Keywords: DNA cation radicals, UV-Vis action spectroscopy, time-dependent DFT calculations, potential energy surfaces

Introduction

Radical reactions play an important role in several biological processes. An area of particular interest has been the reactions triggered in nucleic acids and their components by electrons, photons, and high-energy particles. Nucleobase radicals have been previously generated by electron transfer in the gas phase. We now report on the generation of nucleic acid cation-radicals and their characterization by UV-Vis photodissociation action spectroscopy.

Methods

Electron transfer reduction was used to produce nucleobase cation radicals in a linear ion trap from doubly charged ions formed by electrospray. UV-Vis action spectroscopy was carried out in the 210-70 nm range. Structures were obtained by all-electron molecular dynamics calculations followed by DFT and time-dependent DFT for excitation energies and intensities, including vibronic excitations giving thermal absorption spectra.

Results

Nucleic acid cation radicals of novel type were produced in high yield by electron transfer reduction and dissociation of crown-ether complexes for DNA dimers (AA + 2H), (GG + 2H), (CC + 2H), (GC + 2H), and (CG + 2H). RNA chimeric dimers (A-spacer-A + 2H) were also produced with rigid ring spacers separating the adenosine cation and radical moieties [1]. Action spectroscopy and UV-Vis spectra interpretation revealed that the cation-radicals underwent a conformational collapse associated with hydrogen atom migrations between the nucleobases, forming new isomers. The hydrogen migrations were very facile for adenine and guanine nucleobases, but less so for cytosine. The RNA chimeras with aromatic ring spacers formed zwitterionic structures upon reduction in which the nucleobases maintained the charge and the electron entered the aromatic ring. The formation of zwitterions is consistent with the mechanism and energetics of DNA ionization, as revealed by detailed structure and energy analysis.

Conclusions

Upon electron transfer, single-strand nucleic acids undergo a conformational collapse driven by attractive interactions of nucleobase cations and radicals.

Novel Aspect

Use of UV-VIS action spectroscopy in combination with electron transfer and TD-DFT calculations for structure assignment of transient DNA cation radicals.

References

1. Korn, J. A., Urban, J., Dang, A., Nguyen, H. T. H., Turecek, F. *J. Phys. Chem. Lett.* **8**, 4100 (2017)

MOr.03-4 PHOTOCHEMISTRY OF GLYOXYLATE EMBEDDED IN SODIUM CHLORIDE CLUSTERS, A LABORATORY MODEL FOR TROPOSPHERIC SEA-SALT AEROSOLS

Nina K. Bersenkowitsch - Milan Ončák - Christian van der Linde - Andreas Herburger - Martin K. Beyer

Universität Innsbruck, Ionenphysik und Angewandte Physik, Innsbruck

Keywords: marine aerosols, glyoxylate, photodissociation, $\text{CO}_2^{\bullet-}$

Introduction

As the ocean covers more than 70% of the earth's surface, marine aerosols [1] are highly relevant for atmospheric processes. Photochemical processing of organic matter can be initiated by sunlight, where reactive species can be produced. In the present study, sodium chloride clusters serve as laboratory model for sea salt clusters. Doped with the atmospherically relevant glyoxylic acid [2], the photochemistry is investigated experimentally and supported with theory [3].

Methods

The clusters are produced with an Electrospray Ionization Source, and are then guided into a Fourier Transform Ion Cyclotron Resonance Mass Spectrometer. After mass selection of a specific ion, the clusters get irradiated with a tunable laser system and photodissociation spectra are recorded. Via the intensity of the fragments, photodissociation cross sections are calculated. Theoretical calculations support the results.

Results

Photodissociation spectra from 225-400 nm of $[\text{Na}_n\text{Cl}_{n-2}(\text{C}_2\text{HO}_3)]^+$, $n=5-11$, reveal the fragmentation pathway, where glyoxylate is photolysed into HCO^{\bullet} and a stable carbon dioxide radical anion. This is interesting, as the latter species is known to be metastable in free space, which means that the salt environment stabilizes the molecular anion. Also important is that this happens not only in the deeper UV, but also in the actinic region from 300-400 nm, where sunlight reaches the troposphere. Excited state calculations show that the S_0 and S_1 state are well separated for the isolated glyoxylate anion, which renders internal conversion between these states impossible. As soon as glyoxylate interacts with an ion like H^+ or Na^+ , the potential energy surfaces show a conical intersection, which can be reached with the energy provided from the excitation into S_1 . After funneling into the electronic ground state, the dissociating HCO^{\bullet} may react in the ground state and transfer a hydrogen atom to form formate, which is another observed photochemical fragment.

Conclusions

Photodissociation spectra of sodium iodide clusters of different size doped with a single glyoxylate reveal a fragmentation pathway, where a carbon dioxide anion radical is stabilized in the salt cluster. Quantum chemical calculations show that dissociation along the C–C bond occurs after reaching the S_1/S_0 conical intersection, while this conical intersection is absent in free glyoxylate ions.

Novel Aspect

The photochemistry of glyoxylate embedded in sodium chloride clusters is investigated experimentally, and the results are supported with theoretical calculations.

References

1. Finlayson-Pitts B.J., Pitts J.N., Chemistry of the Upper and Lower Atmosphere. Academic Press, San Diego, 2000.
2. Eugene A.J., Xia S.-S., Guzman M.I., *J. Phys. Chem. A*, 120, 3817 (2016).
3. Bersenkovitsch N.K., Oncak M., van der Linde C., Herburger A., Beyer, M.K., *Phys. Chem. Chem. Phys.*, 20, 8143-8151 (2018).

MO.r.03-5 **PHOTODISSOCIATION SPECTROSCOPY OF M/Z-SELECTED FUNCTIONALISED PROTONATED PYRIDINE IONS: ARE SOME PREDICTIVE PATTERNS EMERGING?**

Benjamin McKinnon ⁽¹⁾ – **Samuel Marlton** ⁽¹⁾ – **Patrick Kelly** ⁽¹⁾ - **Stephen Blanksby** ⁽²⁾ - **Adam Trevitt** ⁽¹⁾

(¹) School of Chemistry, University of Wollongong, Wollongong - (²) Central Analytical Research Facility, Queensland University of Technology, Brisbane

Keywords: laser spectroscopy, ions, ion trap, protonation, photodissociation

Introduction

Protonated N-containing aromatics are ubiquitous fundamental molecules but despite their wide-spread prevalence, there are few predictive models that describe their photostability (and fate). Photostability affects photo-induced proton-transfer, electron-transfer and internal conversion. We are targeting small, protonated N-aromatics to better understand photostability to improve this general understanding.

Methods

A combination of room temperature ion-trap mass spectrometry [1,2], ion-molecule kinetics [3] and tunable UV-Vis pulsed-laser photodissociation is used to measure the photodissociation of *m/z* selected ions formed from electrospray ionisation. Computational results for excited-states and ground state photodissociation pathways assist in assigning the spectra and rationalising the formation of product ions.

Results

The UV photodissociation spectra for ethynyl pyridineH⁺ and formylpyridineH⁺ ions (ortho, meta and para, in each case) are presented for first time. In some cases, the spectra exhibit vibronic detail – these are compared to Franck-Condon simulations – and in all cases multiple photodissociation channels are detected. Ion chemistry kinetics assists with assigning product isomers.

For ethynylpyridiniumH⁺ ions, UV photodissociation leads to either retention of the N atom in the ring (from C₂H₂ elimination) or its expulsion (from HCN elimination). It is not clear what mediates this partitioning. UV photodissociation of formylpyridineH⁺ reveals (at least) three distinct photoproduct channels and these are tracked within the action spectra, which are rather different for each isomer. These different spectra are likely the result of varying excited-state lifetime from to the interplay between “bright” excited states – that dominate the spectra – and “dark” states that dominate the deactivation.

Conclusions

Photodissociation spectroscopy of mass selected ions reveals important new details about their electronic excitation and fate. Although these target molecules are small, there remains many profound unknowns about the fate of their electronic excited states and mechanisms of deactivation. Our study reveals new information about the excited state energy of these species and some plausible mechanisms for their dissociation from UV excitation.

Novel Aspect

New photodissociation spectra are reported for *m/z* selected ethynylpyridiniumH⁺ and formylpyridineH⁺ ions using ion trap MS combined with tunable UV laser irradiation.

References

1. C. S. Hansen, S. J. Blanksby, A. J. Trevitt, Ultraviolet photodissociation action spectroscopy of gas-phase protonated quinoline and isoquinoline cations, *Phys. Chem. Chem. Phys.* **17**, 25882-25890 (2015)
2. C. S. Hansen, S. J. Blanksby, N. Chalyavi, E. J. Bieske, J. R. Reimers and A. J. Trevitt, Ultraviolet photodissociation action spectroscopy of the N-pyridinium cation, *J. Chem. Phys.* **142**, 014301 (2015)
3. C. C. Bright, M. B. Prendergast, P. D. Kelly, J. P. Bezzina, S. J. Blanksby, G. da Silva, and A. J. Trevitt, *Phys. Chem. Chem. Phys.* **19**, 31072-31084 (2017)

MOr.04 - Hyphenated techniques

Chairs: Myeong Hee Moon, Michele Suman

MOr.04-1 **Keynote:** HYPHENATION OF SAMPLING/SAMPLE PREPARATION TO MASS SPECTROMETRY

Janusz Pawliszyn

University of Waterloo, Chemistry, Waterloo

Keywords: Sample Preparation; Ambient MS; SPME; On-site; In-vivo Determinations

Introduction

Since first introduction in 1990, Solid phase microextraction techniques have continually evolved in terms of architecture in order to further push down detection limits while maintaining a clean extraction, portable and easy to handle design. With this portability in mind, the SPME fiber has been appropriately coupled with emerging miniaturized MS technologies for on-site analytical applications. In this way various complex, including living systems can be monitored. The most challenging are biomedical applications. Patient monitoring, metabolomics and biomarkers discovery are an integral part of medical diagnosis and biomedical research. However, tissue analysis remains the bottleneck of such studies due to the invasiveness of presently used sample collection based on biopsy, as well as the laborious and time-consuming sample preparation protocols. In the proposed technology, we have integrated in vivo sampling, sample preparation including global extraction of metabolites together with convenient coupling to analytical instrument resulting in medical biopsy tool.

Methods

This approach based on a small sorbent coated flexible metal fiber format of solid phase microextraction was developed and evaluated during organ treatment and transplantation in rodent or pig model as well as applied to human organs. Other formats will be discussed as well. In this approach, only molecules of interests are extracted onto biocompatible sorbent from the investigated system, however no tissue is not removed from the organ. By using extraction phase with appropriate chemical property and morphology, the analytes are enriched onto the sorbent with minimum of interferences collected. Therefore, the direct coupling of the chemical biopsy probe to mass spectrometry is feasible. Such hyphenation facilitates rapid quantification of the compounds of interests resulting in possibility of close to real time monitoring. This approach has potential in biological and clinical on-site investigations including bedside medical diagnosis.

Results

In the presentation we will discuss diverse SPME-based devices recently developed in our laboratory for the extraction/enrichment of analytes of interest from tissue and small volumes of complex sample matrices, which can be directly coupled with mass spectrometry instruments for rapid analysis. In these applications we use SPME devices containing coatings with matrix compatible morphologies. The approaches tested include Coated Blade Spray (CBS), SPME-transmission mode-direct-analysis-in-real-time (SPME-TM-DART-MS), SPME- microfluidic open interface (SPME-MOI). Total analysis time was typically 5 minutes, but did not exceed 15 minutes and sample volumes ranging between 1-100 μL to the whole organs were used. Sampling/sample-preparation is performed either by spotting the sample onto the SPME-device, or by immersing the SPME-device in a vessel containing the sample or placed directly into tissue. Despite short extraction times, limits of detection in the sub-ng/mL range were obtained, while good accuracy, and linearity were attained for all the studied probes (e.g. therapeutic-drugs, drugs of abuse, and immunosuppressants) in the diverse matrices scrutinized (e.g. urine, plasma,

blood, saliva, brain, muscle and other tissue). In the presentation we will describe exemplary cases in which the mere coupling of SPME to MS is not sufficient to answer relevant analytical questions and the use of strategies that allow for removal of co-extracted interferences or source artifacts, such as Ion Mobility and Multiple Reaction Monitoring with Multistage Fragmentation (MRM3), will also be discussed.

Conclusions

Matrix compatible coatings provide clean extraction suitable for direct coupling to mass spectrometers.

Novel Aspect

SPME biopsy probe hyphenated to MS facilitates on site monitoring even of complex and living systems

References

1. Reyes-Garcés N. et al Anal. Chem. 90, 302 (2018)

MOr.04-2 **SIMULTANEOUS QUANTIFICATION OF SIXTEEN POLYPHENOLS IN 31 LEGUMES BY USING HPLC-MS/MS DYNAMIC-MRM TRIPLE QUADRUPOLE**

Giovanni Caprioli - Massimo Ricciutelli - Sauro Vittori - Astride Franks Kamgang Nzekoue - Gianni Sagratini

School of Pharmacy, University of Camerino, Camerino

Keywords: polyphenols, triple quadrupole, HPLC-MS/MS, legumes

Introduction

Pulses represent one of the most important traditionally dietary component worldwide, supplying proteins, dietary fibers, minerals and vitamins. The Food and Agriculture Organization declared 2016 the International Year of Pulses [1]. The potential health benefits of legumes are attributed to the presence of a wide range of phytochemicals including polyphenols that, according to the literature, are directly associated with antioxidant activity [2].

Methods

Different extraction procedures such as acidic, basic and neutral extraction have been tested at different pH, by combining different temperatures and extraction times and the first one was the best. The separation of the different analytes was performed by using a Kinetex PFP column and the mobile phase was composed by water and methanol both containing 0.1% of formic acid. The analysis has been carried out by using HPLC-MS/MS Dynamic MRM triple quadrupole.

Results

The polyphenols of interest are: shikimic acid, gallic acid, delphinidin, catechin, chlorogenic acid, epicatechin, vanillic acid, caffeic acid, siringic acid, cumaric acid, ferulic acid, 3,5 di-CQA, rutin, kaempferol, cyanidin, quercetin. Chromatographic run is very fast as 16 polyphenols elute within about 10 minutes. The approach used implied Dynamic MRM acquisition mode in which analytes are only monitored while they are eluting from the LC and valuable MS duty cycle is not wasted by monitoring them when they are not expected. The sensitivity of the Dynamic MRM method was higher with respect to the MRM one. The LODs achieved with Dynamic MRM method range from 0.0015 mg/l to 0.015 mg/l and LOQs ranged from 0.001 mg/l to 0.05 mg/l. Considering linearity, the obtained coefficients of correlation were in the range 0.9907-0.9998 for the MRM method and 0.9934-1 for the Dynamic MRM method.

Conclusions

The application to real samples indicated that black beans (458.64 mg/kg) and ruviotto beans (189.02 mg/kg) displayed the highest amount of total polyphenols. This study highlights that the content of polyphenols is related with the dark color of coat pulses and that certain polyphenols are characteristic in specific pulses variety. These aspects may be useful for increasing the health knowledge of legumes and for proposing new functional foods.

Novel Aspect

A new analytical method that uses HPLC-MS/MS and triple quadrupole in Dynamic MRM modewas developed for the analysis of sixteen polyphenols in pulses.

References:

1. Food and Agriculture Organization (FAO), 2016. <http://www.fao.org/pulses-2016/en/>
2. Adebamowo C.A., Cho E., Sampson L., Katan M.B., Spiegelman D., Willett W.C., Holmes M.C. *International Journal of Cancer* 114, 628-633 (2005)

MOr.04-3 ON-LINE COUPLING OF SURFACE PLASMON RESONANCE TO AMBIENT MASS SPECTROMETRY

Huwei Liu - Yiding Zhang - Xianjiang Li - Honggang Nie - Yu Bai

Peking University, College of Chemistry and Molecular Engineering, Beijing

Keywords: Ambient mass spectrometry, surface plasmon resonance, hyphenation method, molecular interaction

Introduction

Surface plasmon resonance (SPR) is a powerful technique widely utilized in the studies based on interaction processes such as cellular analysis [1] and clinical researches [2]. Despite its versatility, there still exist some intrinsic deficiencies due to the principle of SPR, among which the insufficiency in qualitative analysis remains a significant one. Generally, the information acquired by conventional SPR instruments is inadequate for the identification of the samples, constraining the application of SPR in the researches dealing with unknown species such as ligand fishing [3] or active ingredient identification [4]. Consequently, the combination of SPR with other analytical approaches that can provide molecular information has been highly desired [5-9].

In this work, a newly developed interface to couple SPR with MS using spray tip and DART was reported. The SPR flow solution was led into a sprayer, whose tip was put in middle of the DART outlet and the MS inlet, thus the solution could be sprayed and ionized. Proof-of-concept experiments using 4 chemicals dissolved in 3 types of salt-containing solvents were conducted, showing the practicability of this interface by observing the consistent change of the signals of MS and SPR. The performance of this interface and its salt tolerance suggested its high potential in the screening of the unknown interaction under physiological conditions using SPR-MS.

Methods

SPR sensor chip modification. The sensor chip was modified with HSA using common EDC/NHS coupling method. The sensor chip was first cleaned by immersing it in a solution blended by ammonia solution, hydrogen peroxide and water with the volume ratio of 1:1:5 and being kept at 90 °C for 10 min. After installing the cleaned chip into the SPR analyzer, the following modification process was carried out by successively injecting and flowing: a. 3-MPA with TCEP, for 60 min; b. EDC and NHS for 20 min; c. HSA for 12 h; d. ethanolamine for 10 min. The change of SPR angle before and after the HSA modification proved that the protein had been successfully immobilized onto the sensor chip.

SPR-AMS experiments. Acetaminophen, metronidazole, quinine and hippuric acid were selected to be the model analytes, with PBS, TBS and EBSS chosen as the representative solvents to simulate the physiological conditions that might be used in usual SPR experiments. There were all together 12 kinds of sample solutions. A single run was started by commencing SPR and MS data acquisition simultaneously. For the first 2 minutes, pure solvent was delivered. At the beginning of the 3rd minute, the sample solution started to flow, and lasted for 4 minutes. During these 4 minutes, the analyte might bind to the HSA immobilized on the sensor surface, leading to a change in the SPR signal. From the 6th minute, again the pure solvent was delivered, and the analyte would be washed away. After the SPR signal fell back to its baseline and the MS signal became stable, the data acquisition was stopped. The raw SPR data was denoised by fast Fourier transformation smoothing before further analysis.

The solution flow rate of all of the runs was 0.20 mL/min. The nebulizing gas was nitrogen at a pressure of 15 psi. The SPR analyzer was set at constant angle mode recording the reflected light intensity at the SPR angle. For the DART-MS ion source, the temperature was set at 450°C with helium as working gas. The MS parameters were: capillary voltage, off; dry gas temperature, 325 °C; dry gas flowrate, 5 L/min. The range of m/z between 100 and 400 was scanned with a maximum ion accumulation time of 200 ms. MS was operated at positive mode when analyzing solutions containing acetaminophen, metronidazole or quinine, and negative mode for solutions of hippuric acid.

Results

Optimization of the DART conditions. As preliminary experiments, the following parameters that might affect the DART ionization and the MS analysis were optimized: DART source temperature, the vertical distance between the spray tip and the DART-outlet-MS-inlet plane (d_3), the total horizontal distance between the DART outlet and the MS inlet (d_1), the flow rate of the solution and the nebulizing gas pressure for the spray. As these parameters were independent, a series of experiments based on the orthogonal experimental design were conducted, using 25 mmol/L acetaminophen dissolved in PBS as the sample.

For d_3 , the temperature and the flow rate, higher values would lead to better signals, and for the nebulizing gas pressure, a lower value would be better. The values of d_1 showed no clear tendency. These effects could be explained by the principle of the interface. As the interface was spray-based, a better spray condition was of great concern. The droplets in the spray would get smaller and better distributed as they flew away from the spray tip, so a larger d_3 would be beneficial. Also, as the analytes were dissolved in the solutions and the sprayed droplets, the ionization would consist of the evaporation and desolvation process, thus a higher source temperature would be more productive. For the nebulizing gas pressure, if it was too high, the droplets would be blown away too fast from the DART ionizing stream, so a lower pressure would be better. The flow rate meant the amount of sample, and a higher value led to a stronger signal. The value of d_1 showed least importance, thus it was selected considering other parameters.

The horizontal distance between the spray tip and the MS inlet (d_2) could also affect the results. However, as there was a restrictive requirement that d_2 must be less than d_1 , the value of d_2 could not be included in the orthogonal experiments and it was optimized separately. Generally, a smaller d_2 could allow more ions to go into the MS inlet, but it might also lead to the sample residue at the MS inlet. Also, a larger d_3 would result in more residues. With consideration to the signal intensity and the sample residue contamination, the following parameters were selected: DART source temperature, 450 °C; d_1 , 2.0 cm; d_2 , 10 mm, d_3 , 6 mm; flow rate, 0.20 mL/min; nebulizing gas pressure, 15 psi. Additionally, helium was used as working gas for all of the sample solutions.

Performance of the SPR-DART-MS experiments.

As mentioned in the experimental section, there were 12 sample solutions. It was undoubtedly proven that the as-designed spray and DART based interface could be used successfully and

practically for on-line coupling of SPR with MS, as the MS signal changed consistently with the SPR signal. For the first 2 minutes during which the pure solvent flowed through the SPR analyzer, SPR signal stayed at its baseline, so did the MS signal. As the solution containing analytes started to flow at the 3rd minute, the SPR signal started to increase as the analyte bound with the HSA immobilized on the sensor chip. The MS signal also increased, as the excess analyte in the solution got sprayed and ionized. From the 6th minute, pure solvent flowed again, and the SPR signal changed to its original values because of the dissociation of the bound analytes, and the MS signal also decreased.

There is an observable lag of approximately 1 min between the increase of the MS signal and the injection of the sample solution. This was mainly because of the volume of the SPR flow cell and the connecting tubes. All of these volumes between the pump and the MS would act like the dead volume similar in a chromatography process, delaying the change of the MS signal. Besides the dead volume, there were another two reasons explaining the delay of this decrease. First, there remained analytes bound on the sensor chip, the elution of these analytes would last for some time, thus the MS signal would remain high during all this elution process. Second, as the aqueous solutions containing non-volatile salts were used, the evaporation of the sprayed droplets would often be incomplete, resulting in the residue of the analytes in the vicinity of the spray tip and the MS inlet, keeping the MS signal from decreasing. Notwithstanding this residue, it demonstrated the necessity and advantage of the DART-based interface, that it could prevent the residue from contaminating conventional ESI source.

Salt tolerance of the DART interface.

For most analytes and solvents, the interface showed a high tolerance of the non-volatile salts. For the 3 analytes in positive mode, both of PBS and EBSS showed even higher signal intensity than water, indicating that the low-molecular-weight salts might help the proton transfer to form positive species. TBS showed lower signal intensity, which could be explained by the competitive ionization between the analyte and tris contained in TBS, which was supported by the strong tris peak ($m/z = 122$, $[M+H]^+$) observed in most of the MS spectra using TBS. For the analyte in negative mode, all of the salt containing solvents showed a depression in signal intensity, suggesting the salts' negative effects in negative mode. These results proved the as-designed interface had an adequate tolerance of the non-volatile salts, and should suffice the applications for direct on-line coupling of SPR with MS using salt-containing solutions.

Conclusions

A novel interface relying on DART ionization and spray tip for direct and on-line coupling of SPR with MS was designed and developed. Experimental results using different model compounds and buffers exhibited the feasibility of the online combination of SPR and MS. This system showed a high tolerance of the non-volatile salts during the test, suggesting its potential utilizations in SPR-MS experiments, especially when physiological conditions should be needed. This interface provided a new approach for online coupling of SPR-MS, and would contribute to the fast detection of the interaction between substrate and the ligand and simultaneous identification of the ligand [10].

Novel Aspect

A novel interface for direct and on-line coupling of SPR to AMS was developed. A spray tip connected with the outlet of the SPR flow solution was conducted as the sampling part of the DART-MS, with which the online coupling interface of SPR-MS was realized. Four model samples, acetaminophen, metronidazole, quinine and hippuric acid dissolved in three kinds of common buffers were used in the SPR-DART-MS experiments for performance evaluation of the interface and the optimization of DART conditions. The results showed consistent signal changes and high tolerance of non-volatile salts of this SPR-MS system, demonstrating the feasibility of the interface for on-line coupling of SPR with MS and the potential application in the characterization of interaction under physiological conditions.

References

1. Abadian P. N., Kelley C. P., Goluch E. D., Anal. Chem., 86, 2799-2812 (2014).
2. Mariani S., Minunni M., Anal. Bioanal. Chem., 406, 2303-2323 (2014).
3. Peng M., Zhang Y., Shi S., Peng S. J., Chromatogr. B., 940, 86-93 (2013).
4. Zhang Y., Shi S., Guo J., You Q., Feng D. J., Chromatogr. A, 1293, 92-99 (2013).
5. Krone J. R., Nelson R. W., Dogruel D., Williams P., Granzow R., Anal. Biochem., 244, 124-132 (1997).
6. Nelson, R. W.; Krone, J. R. J. Mol. Recognit., 12, 77-93 (1997).
7. Nedelkov, D.; Nelson, R. W. Trends. Biotechnol., 21, 301-305 (2003).
8. Visser, N. F. C.; Heck, A. J. R. Expert Rev. Proteomics, 5, 425-433 (2008).
9. Stigter, E. C. A.; de Jong, G. J.; van Bennekom, W. P. TrAC-Trend. Anal. Chem., 45, 107-120 (2013).
10. Zhang, Y.; Li X.; Nie, H., Yang, L.; Li, Z.; Bai, Y.; Niu, L.; Song, D.; Liu, H., Anal. Chem., 87, 6505-6509 (2015).

MOr.04-4 **ONLINE HIGH PRESSURE PROTEOLYTIC BIOSENSOR- MASS SPECTROMETRY (PROTEX-SPRMS): A NEW POWERFUL TOOL FOR PROTEIN EPITOPE AND AFFINITY IDENTIFICATION AND ANTIBODY CHARACTERIZATION**

Loredana Lupu⁽¹⁾ - **Pascal Wiegand**⁽¹⁾ - **Yannick Baschung**⁽¹⁾ - **Friedemann Völklein**⁽²⁾ - **Alexander Lazarev**⁽³⁾ - **Katja Schmitz**⁽⁴⁾ - **Michael Przybylski**^(1,2)

⁽¹⁾ Steinbeis Center for Biopolymer Analysis & Biomedical Mass Spectrometry, Marktstraße 29, 65428 Rüsselsheim am Main, Germany - ⁽²⁾ Rhein Main University, Institute of Microtechnologies, 65428 Rüsselsheim am Main, Germany - ⁽³⁾ Pressure Biosciences Inc., Boston - ⁽⁴⁾ Techn. University Darmstadt, Department of Organic Chemistry and Biochemistry, 64287 Darmstadt, Germany

Keywords: Online-Biosensor-MS combination; T-SPR biosensor; continuous flow microfluidics; antibody epitope identification; epitope peptide affinity determination.

Introduction

Analysis of biomolecular interaction epitopes is a key step in the development of therapeutic antibodies, biomedical peptide and protein biomarkers and molecular vaccines. Bioaffinity analysis using biosensors has become an established technique for detection and quantification of biomolecular interactions. However, a principal limitation of biosensors is their lack of providing chemical structure information of affinity-bound ligands. Proteolytic excision/extraction (PROTEX-MS) is a major techniques for mass spectrometry based elucidation of protein- ligand interactions, but affinity-MS analysis does not provide quantitative affinity data. A new online proteolytic excision SPRMS combination (PROTEX-MS) enables simultaneous structure identification and affinity quantification of antibody epitopes from a protein- ligand complex immobilized on a gold chip.

Methods

Key tools of the online proteolytic excision SPRMS Epitope Analyzer are (i), a new SPR biosensor that enables the application of microfluidic systems suitable for direct ESI-MS that is not sensitive to buffer change normally required for SPR (T-SPR); (ii), an integrated, automated interface that provides elution and sample concentration for direct MS analysis of the eluted epitope peptides [1]; and (iii), a new high pressure proteolytic system, high pressure barocycler digestion. The high pressure digestion enables substantially enhanced speed, efficiency and automation of protease digestion and proteolytic affinity-extraction, compared to conventional digestion.

Results

Initial validation of the new PROTEX-To validate SPR-MS system, analysis on a well-defined affinity pair has been studied. Results from the tests on anti-Beta Amyloid antibody interacting with Abeta 1-40 and 1-16 protein provided us with a K_D values of 0.628 μM and 0.796 μM and with the protein/peptide structural information. The broad application potential of the PROTEX-SPR-MS Epitope Analyzer has been ascertained by identification of protein epitope structures from therapeutic and pathophysiological antibodies, such as an unusual mixed-disulfide antibody epitope of the rheumatic target protein, HLA-B27; and the interaction site identification of chaperone complexes of lysosomal enzymes [2, 3]. Interaction epitopes as diverse as antigen-antibody and lectin- carbohydrate complexes [4, 5], with affinity binding constants (K_D) from milli- to nanomolar ranges can be directly analyzed. Applications amenable with the PROTEX-SPR-MS Epitope Analyzer include affinity-based biomarker identification; identification of protein and peptide epitopes; precise antibody affinity determinations; direct label-free antigen quantification; and clinical epitope analyses of pathophysiological antibodies, suitable for clinical apheresis and hyposensitization of patients after immunogenic protein therapy [4-6].

Epitopes from antisera of Fabry disease (FD; α -Galactosidase deficiency) patients were identified using the PROTEX-SPR-MS combination. Tryptic peptide mixtures were loaded onto the SPR-MS; after washing out nonbinding peptides, epitope peptides were eluted with 0.01% TFA into the affinity interface, ESIMS was performed with a Waters QuattroUltima MS. The epitope peptide, α -Gal(309-332) was identified from 2 FD patients and revealed high affinity (K_D , 39 nM) comparable to the full length enzyme (K_D , 16 nM). Linear, stabilized epitope peptides were prepared by SPPS.

Conclusions

A number of applications are amenable with the new online PROTEX-SPR-MS epitope analyzer, including affinity-based biomarker identification, identification of protein and peptide epitopes, precise antibody affinity determinations, and direct label-free antigen quantification. In addition, the clinical application potential was shown in epitope elucidation from therapeutic and pathophysiological antibodies of lysosomal enzymes, with affinity constants (K_D) amenable from milli- to nanomolar ranges.

Novel Aspect

A new powerful tool for epitope and affinity analysis of biomolecular interactions was developed. The PROTEX- SPRMS Epitope Analyzer provides a breakthrough (i) for proteolytic digestion of proteins; and (ii), direct affinity determination in combination with MS.

References

1. Przybylski M, Lupu L, Voelklein F (2018) Eur. Pat. Appl.
2. Iurascu M, Belaunzar OM, Cozma C, Petrausch U, Renner C, Przybylski M (2016) J. Am. Soc. Mass Spectrom. 27, 1105-1112.
3. Moise A, Maeser S, Rawer S, Eggers F, Murphy M, Bornheim J, Przybylski M (2016) J. Am. Soc. Mass Spectrom. 27, 1071-1078.
4. Paraschiv G, Vincke C, Czapplewska P, Manea M, Muyldermans S, Przybylski M (2013) J. Mol. Recogn. 26, 1-9
5. Kukacka Z, Iurascu M, Lupu L, Rusche H, Murphy M, Altamore L, Borri F, Maeser S, Papini AM, Hennermann J, Przybylski M (2018) ChemMedChem., DOI 101002 201800094.

MOr.04-5 **CHEMICAL PROFILING OF THE HUMAN SKIN SURFACE FOR MALARIA VECTOR CONTROL USING A NON-INVASIVE SOLVENT MODIFIED SORPTIVE SAMPLER WITH GCxGC-TOFMS AND UPLC-IMS-HRMS**

Madelien Wooding⁽¹⁾ - **Tyren Dodgen**⁽²⁾ - **Marc Bouwer**⁽³⁾ - **Yvette Naudé**⁽¹⁾ - **Egmont Rohwer**⁽¹⁾

⁽¹⁾ University of Pretoria, Department of Chemistry, Pretoria - ⁽²⁾ Microsep, Application Specialist, Johannesburg - ⁽³⁾ University of Pretoria, Forestry and Agricultural Biotechnology Institute, Pretoria

Keywords: GC×GC-MS, UPLC-IMS-HRMS, human skin compounds, sorptive sampling, malaria vector control

Introduction

Odour mosquito lures are currently being used as part of integrated vector control strategies in the fight against malaria. Variation in inter-human attractiveness to mosquitoes, as well as the preference of mosquitoes to bite certain regions on the human host are possible avenues for identifying lead compounds as potential attractants[1-4]. These may prove more effective in lure and kill strategies.

Methods

Two different areas of the human skin, namely wrist and ankles, were compared following a 1 hour sampling period. Skin surface compounds were concentrated using a solvent modified polydimethylsiloxane (PDMS) sampler. Sampling was followed by analyses with GC×GC-TOFMS and UPLC-IMS-HRMS. The sampler was thermally desorbed directly in the inlet liner of a GC system, while for LC analyses the PDMS sampler was solvent back-extracted prior to analyses.

Results

The skin profiles for the ankle and wrist groups were compared using principal component analyses. Compounds responsible for the differences in the chemical profiles were tentatively identified based on a comparison of sample mass spectra to that of the NIST14 library for the GC analyses. Accurate mass, isotope fit values and fragmentation patterns were compared to online ChemSpider databases for the LC analyses. Binary compare using an OPLS-DA Model highlighted the differences between the ankle and wrists groups for the LC-MS data sets. Viburtinaland Menaquinol, as well as Allodeoxycholic acid and *N,N*-Diethylbenzeneacetamide were tentatively identified, using the Human Metabolome Database. These form plausible lead markers for the wrist and ankle groups, respectively. A comprehensive list of collision cross section data of human skin chemicals was compiled for LC amenable markers.

Conclusions

In this study, methods were developed to chemically profile the human skin surface with the aim of determining the compounds responsible for attracting the mosquito vector to its human host. A PDMS sorptive sampler was modified with solvent to give improved extraction of polar analytes from the skin surface. This multi-approach gave a comprehensive human skin chemical profile consisting of volatile, semi-volatile, and non-volatile compounds.

Novel Aspect

Skin surface compounds were concentrated using a solvent modified, non-invasive, in-house developed PDMS sampler for GC-MS and LC-MS. CCS library for skin surface compounds.

References

1. Takken, W. and B.G. Knols, *Annual review of entomology*, 44(1), p. 131-157 (1999).
2. Hallem, E.A., et al., *Nature*, 427(6971), p. 212-213 (2004).
3. Verhulst, N.O., et al., *PloS one*, 6(12), p. e28991(2011).
4. Clements, A.N., *The biology of mosquitoes: sensory reception and behaviour*, Cambridge, UK: Cambridge University Press(1999).

MOr.05 - Volatile molecules & nose-space analysis

Chairs: Egmont Rohwer, Giorgio Mellerio

Monday

MOr.05-1 **Keynote: EXPLORING FOOD AROMA RELEASE AND PERCEPTION USING DIRECT-INJECTION MASS SPECTROMETRY**

Jonathan Beauchamp

Fraunhofer Institute for Process Engineering and Packaging IVV, Department of Sensory Analytics, Freising, Germany

Keywords: DIMS, aroma, nosespace, olfaction, perception

Introduction

Direct-injection mass spectrometry (DIMS) –in the forms of proton-transfer-reaction MS (PTR-MS), atmospheric pressure chemical ionization MS (APCI-MS), and selected ion flow tube MS (SIFT-MS) – has been increasingly implemented over the past two decades to explore food flavor release and perception. This keynote lecture will review novel food flavor insights that have been achieved using on-line DIMS techniques, focusing on applications in *in vivo* analysis.

Methods

Food flavor perception is a complex process that depends on manifold parameters, including the physicochemical properties of the flavor compound, the food matrix, saliva, oral and nasal anatomies, as well as the mucosa lining the airways, amongst other factors [1]. The real-time analytical capability of DIMS techniques has opened up new opportunities for characterizing the fast dynamics of flavor release and relating this to perceived sensory impressions.

Results

In vivo analyses using DIMS was pioneered around two decades ago by Taylor and Linforth using APCI-MS [2,3,4] and offers an opportunity to explore parameters affecting flavor release and perception systematically. Analyses involve sampling exhaled breath via the nose – termed nosespace analysis – to capture volatile aroma compounds that are released from the food (or beverage) matrix during oral processing and enter the nasal cavity during swallowing, eliciting perceived flavor impressions. This approach has been used, for example, to assess aroma perception in relation to taste attributes [5] and textural and physiological parameters [6], or olfactory acuity in general [7]. DIMS has also been applied to generate sensory profiles of food products based on non-destructive headspace analysis, as a means to characterize a product's sensory attributes without the (explicit) need of a trained human sensory panel [8]. New developments and peripheral equipment have extended the range and versatility of DIMS systems in food flavor research [9,10].

Conclusions

Nosespace analyses via DIMS techniques has provided many insights into flavor perception that would not have been achievable using conventional mass spectrometry in combination with classical gas chromatography (GC-MS). The shortcomings of DIMS techniques tend to be their lack of unequivocal identification, but this requirement is obviated when using a targeted approach to monitor fast processes of individual flavor compounds that are known *a priori*.

Novel Aspect

On-line direct-injection mass spectrometry is a versatile and powerful tool for investigating the dynamic processes involved in flavor release and perception.

References

1. Buettner A., Beauchamp J., *Food Qual. Pref.*, 21, 915-924 (2010).
2. Baek I., Linforth R. S., Blake A., Taylor A. J., *Chem. Senses*, 24, 155-160 (1999).
3. Linforth R., Taylor A. J., *Journal of Agricultural and Food Chem.*, 48, 5419-5423 (2000).
4. Taylor A. J., Linforth R. S. T., Harvey B. A., Blake A., *Food Chem.*, 71, 327-338 (2000).
5. Davidson J. M., Linforth R. S. T., Hollowood T. A., Taylor A. J., *J. Agric. Food Chem.* 47, 4336-4340 (2000).
6. Mestres M., Kieffer R., Buettner A., *J. Agric. Food Chem.*, 54, 1814-1821 (2006).
7. Beauchamp J., Scheibe M., Hummel T., Buettner A., *Chem. Biodivers.*, 11, 619-638 (2014).
8. Lindinger C., Labbe D., Pollien P., Rytz A., Juillerat M. A., Yeretian C., Blank I., *Anal. Chem.* 80, 1574-81 (2008).
9. Beauchamp J., Herbig J., *The Chemical Sensory Informatics of Food: Measurement, Analysis, Integration*, ed B Guthrie, et al.: American Chemical Society) pp 235-51 (2015).
10. Le Quéré J.-L., Gierczynski I., Sémon E., *J. Mass Spectrom.* 49, 918-28 (2014).

MOr.05-2 **HIGH-PRESSURE PHOTOIONIZATION TIME OF FLIGHT MASS SPECTROMETRY COMBINED WITH DYNAMIC PURGE-INJECTION FOR RAPID ANALYSIS OF VOLATILE METABOLITES IN URINE**

Hua Lei⁽¹⁾ - **Lei Hua**⁽¹⁾ - **Enyou Li**⁽²⁾ - **Haiyang Li**⁽¹⁾

⁽¹⁾ Dalian Institute of Chemical Physics, Chinese Academy of Sciences, CAS Key Laboratory of Separation Sciences for Analytical Chemistry, Dalian - ⁽²⁾ The First Affiliated Hospital of Harbin Medical University, Department of Anesthesiology, Haerbin

Keywords: high-pressure photoionization mass spectrometry, dynamic purge-injection, urine analysis, volatile organic metabolites

Introduction

Small molecule metabolites are widely used as biomarkers in the research field of metabolomics for disease diagnosis and exposure assessment. As a readily available biofluid containing plenty of volatile organic metabolites (VOMs), urine is ideal for non-invasive metabolomics analysis.[1] However, there is still lack of rapid analysis method for VOMs in urine.

Methods

A dynamic purge-injection apparatus was designed by coupling of a thermostated water bath and a bubbling bottle. The urines were firstly acidized or alkalized by simply adding a small amount of H₂SO₄ or NaOH. The VOMs were released into the gaseous phase by bubbles bursting at the gas-liquid interface, which were taken into the high-pressure photoionization (HPPI) mass spectrometer for analysis.[2]

Results

Various types of metabolites, such as ketones, alcohols, acids, sulfides, pyrroles and amines were detected directly by simple acidification or alkalization of urines. Nitrogen-containing compounds, especially polar amines, could be ultra-sensitively measured without any derivatization. Interestingly, a new biomarker 2,5-dimethylpyrrole was exclusively found in the smoker's urine sample besides toluene. Five valuable metabolites, trimethylamine, toluene, styrene, p-xylene, and 2,5-dimethylpyrrole, were well quantified with the limits of detection down to 0.02–0.1 µg L⁻¹, and the relative standard deviation (R.S.D.) ranged from 0.42% to 14.4%. Being highly sensitive and humidity-friendly, the whole analytical procedure was easily operated in less than 6 min.

Conclusions

The coupling of HPPI-TOFMS with dynamic purge-injection was demonstrated as a novel tool for

rapid nontarget disease biomarkers screening or target monitoring of specific compounds through the investigation of volatile metabolites in urine. The highly sensitive and humidity-friendly system might aid in clinical diagnosis for development of diagnostic tools for disease identification.

Novel Aspect

A powerful and sensitive direct mass spectrometric technique for high-throughput analysis of volatile organic metabolites in urine.

References

1. Calejo I., Moreira N., Araujo A.M., Carvalho M., Bastos Mde L., de Pinho P.G., *Talanta*, 148, 486-493 (2016).
2. Wang Y., Hua L., Jiang J., Xie Y., Hou K., Li Q., Wu C., Li H., *Analytica Chimica Acta*, 1008, 74-81 (2018).

MOr.05-3 NOSE-SPACE OF DARK CHOCOLATES USING PTR-TOF-MS AND LINK TO FLAVOUR PERCEPTION THROUGH SIMULTANEOUS TEMPORAL DOMINANCE OF SENSATIONS (TDS)

Zoé Deuscher⁽¹⁾ - Isabelle Andriot⁽¹⁾ - Sylvie Cordelle⁽²⁾ - Marie Repoux⁽³⁾ - Renaud Boulanger⁽⁴⁾ - H  l  ne Labour  ⁽⁵⁾ - Pascal Schlich⁽¹⁾ - Jean-Luc Le Qu  r  ⁽⁶⁾

⁽¹⁾ INRA, CSGA, Dijon - ⁽²⁾ CNRS, CSGA, Dijon - ⁽³⁾ Valrhona, Valrhona, Tain-l'Hermitage - ⁽⁴⁾ CIRAD, Qualisud, Montpellier - ⁽⁵⁾ AgroSup Dijon, CSGA, Dijon - ⁽⁶⁾ INRA, Centre des Sciences du Go  t et de l'Alimentation (CSGA), Dijon

Keywords: chocolate – nosespace – PTR-ToF-MS – Temporal Dominance of Sensations (TDS)

Introduction

Aroma of dark chocolate depends on process and cocoa origin and variety. A sensory analysis of 206 standardized chocolates produced from various cocoa beans classified them in four sensory categories. These categories were confirmed in a PTR-MS analysis of the chocolates volatilome [1]. The objective here was to study the nosespace of a subset of chocolates simultaneously with their temporal profile to better explain the sensory categorization at a perception level.

Methods

A Temporal Dominance of Sensations (TDS) evaluation of 8 chocolates (2 selected per category) was done in triplicate by 12 subjects while the aroma released in their nose were simultaneously collected and injected into a PTR-ToF-MS. The two sets of data were analyzed conjointly by defining an index of abundance of each detected aroma compound while a given attribute was dominant: the Abundance While Dominance (AWD) index [2].

Results

TDS is able to dynamically capture multidimensionality of perception. The obtained TDS curves clearly differentiated the chocolates that were regrouped by sensory category as revealed by a principal component analysis (PCA). The dynamics of perception have been considered only scarcely in comparison to simultaneous in vivo aroma release over time measured by on-line mass spectrometry (nosespace). Computation of the AWD indices at individual level allowed to statistically assess the differences between the products over subjects and replicates and assessed statistically the relationships between the two sets of data. Through correspondence analyses (CA) some relationships between certain aroma compounds and the sensory attributes expected to be related to them were found. Although in previous studies pairing nosespace and TDS various temporal links could be proposed [3], no clear relationships could be safely

established due to the fact that the conclusions were mainly based on a descriptive analysis of the data conducted at panel level.

Conclusions

Pairing nosespace with Temporal Dominance of Sensations evaluation of dark chocolates categorized in four sensory groups provides meaningful data that can be analyzed at individual level thanks to the AWD indices. Descriptive multivariate analyses of these AWD indices gave interesting clues on the relationships between the aroma compounds released in mouth and their expected perceived sensory attributes.

Novel Aspect

Pairing nosespace PTR-ToF-MS with Temporal Dominance of Sensations allows assessing the relationships between aroma compounds released in vivo and their perception.

References

1. Deuscher Z., Andriot I., Sémon E., ..., Le Quéré J.L., *SMMAP2017*, Dineyland Paris, October 2-5 (2017).
2. Schlich P., Thomas A., Visalli M., Labarre D., Sémon E., Le Quéré J.L., in *Flavour Science: Proceedings of the XIV Weurman Flavour Research Symposium*, (Taylor A.J., Mottram D.S., Eds.), Context Products Ltd., Leicestershire (UK), 327-332 (2015).
3. Déléris I., Saint-Eve A., Dakowski F., Sémon E., Le Quéré J.L., Guillemin H., Souchon I., *Food Chem.*, 127(4), 1615-1624 (2011).

MOr.05-4 **IN NOSE-SPACE ANALYSIS THROUGH PROTON TRANSFER REACTION MASS SPECTROMETRY: INVESTIGATING VARIABILITY IN FLAVOR RELEASE AND PERCEPTION DURING CHEWING GUM CONSUMPTION**

Michele Pedrotti⁽¹⁾ - **Iuliia Khomenko**⁽²⁾ - **Vincenzo Fogliano**⁽³⁾ - **Franco Biasioli**⁽²⁾

⁽¹⁾ Foundation Edmund Mach, Wageningen University FQD, Wageningen - ⁽²⁾ Foundation Edmund Mach, DQAN, San Michele all'Adige - ⁽³⁾ Wageningen University, FQD, Wageningen

Keywords: nose-space analysis, chewing gum, in vivo flavor analysis, PTR-MS

Introduction

Proton Transfer Reaction Mass Spectrometry has been described as particularly suited to provide “in vivo” on-line flavor monitoring during food consumption [1]. In the last years, the technique has been applied to different food matrixes [2-4]. Due to analysis rapidity and high sensitivity, the technique has become a reference tool for the investigation of the complex phenomena related to flavor perception [5].

In this work we describe a PTR-MS based approach to measure in-vivo nose-space VOCs concentration during chewing gum consumption coupled to discontinuous time intensity sensory evaluation.

Methods

Two different studies were conducted on chewing gums to monitor respectively differences in aroma release due to two different gum ingredients formulations and the effect of physiological parameters (oral cavity volume, salivary flow and papillary count), gender and ethnicity in flavor perception. Sensory evaluation on flavor intensity and sweetness was performed simultaneously with nose-space analysis with two different PTR-MS machines (PTR-ToF-MS and PTR-QiTOF-MS, Ionicon Analytik, Austria). A total of 40 panelists participated in the two studies.

Results

In the first study, a significant difference in total aroma release (p.value < 0.05 ANOVA) was

observed between the different chewing gum formulations. In particular, it was observed that these differences were compound dependent and that PTR-MS signals are associated to a high variability reflecting panelists' physiological differences.

In the second study, ethnicity was found to have a significant effect on both in nose-space concentration and sensory perception. For different mass peaks associated to mint flavor compounds, Chinese panelists exhibited higher levels than European ones (p.value < 0.05 Welch's t-test) generally after 90 seconds of consumption and after the gum was removed from the mouth. The same trend was found in sensory perception both for flavor and sweetness attributes. No differences were found between panelists of different gender in both flavor release and perception.

Conclusions

Real time mass spectrometry analysis by PTR-MS can detect and monitor in-vivo volatile organic compounds release during food consumption.

The results of the studies suggest a multimodal effect of aroma release on both flavour and sweetness intensity. Individual differences due to physiological, biochemical and physiocochemical phenomena may have a relevant effect in aroma perception and thus in sensory flavor perception. Although, it cannot be excluded a cultural bias due to an ethnocentric experimental environment, result suggest that higher in-nose concentration in Chinese judges may explain their higher flavor and sweetness perception.

Novel Aspect

Real time direct injection mass spectrometric analysis coupled with dynamic sensory methods is a powerful tool to investigate flavor perception mechanisms and highlight physiological and cultural biases.

References

1. Biasioli F., Gasperi F., Yeretian C., Märk T.D., *Trends in Analytical Chemistry*, , 30, 968–977 (2011).
2. Romano A., Cappellin L., Ting V., Aprea E., Navarini L., Gasperi F., et al., *International Journal of Mass Spectrometry*, 365–366, 20–27 (2014).
3. Heenan S., Soukoulis C., Silcock P., Fabris A., Aprea E., Cappellin L., Märk T.D., Gasperi F., Biasioli F., *Food Chemistry*, 131, 477 (2012).
4. Fiches G., Saint-Eve A., Jourden S., Délérís I., Brunerie P., Souchon I., *Flavour and Fragrance Journal*, 31, 31-40 (2016).
5. Avison J., *Journal of Agricultural and Food Chemistry*, 61, 2070-2076 (2013).

MOr.05-5 **METHODOLOGY FOR QUALITATIVE AND QUANTITATIVE ANALYSIS OF VOLATILE COMPOUNDS FROM COMPOSITE MATERIALS AT ELEVATED TEMPERATURES**

Célia Lourenço⁽¹⁾ - **Sarah Bergin**⁽¹⁾ - **Daniel Francis**⁽¹⁾ - **Stephen Staines**⁽¹⁾ - **Jane Hodgkinson**⁽¹⁾ - **Chris Walton**⁽²⁾ - **John Saffell**⁽³⁾ - **Ralph P. Tatam**⁽¹⁾

⁽¹⁾ Cranfield University, Cranfield University/Engineering Photonics, Cranfield, Bedfordshire - ⁽²⁾ Cranfield University, Cranfield University/Centre for Environmental and Agricultural Informatics, Cranfield, Bedfordshire - ⁽³⁾ Alphasense, Alphasense, Essex

Keywords: Volatile organic compounds, composites, gas sensors, TD-GC-MS

Introduction

The increasing use of composite materials on the aircraft structure led to the need to investigate their performance under on-board operating conditions and non-standard conditions such as the

onset of fire [1]. The investigation of possible effects on the on-board air quality arising from the introduction of new composite materials on the aircraft structure, particularly at elevated temperatures where Volatile Organic Compounds (VOCs) may be released is evaluated here.

Methods

The described methodology allowed the characterisation of composite materials inside an air-tight tube furnace operated up to 250°C. The real-time detection of released gases was accomplished using commercial off-the-shelf (COTS) gas sensors, in addition, stainless steel thermal desorption (TD) tubes were used and analysed by gas chromatography-mass spectrometry (GC-MS) combined with thermal desorption (TD-GC-MS) for both qualitative and quantitative analysis.

Results

Real-time detection of released gases (CO, NO, SO₂, NO₂, CO₂, O₂) and volatiles was performed over a temperature ramp from 22°C with dwell points at 70°C, 150°C, 200°C and 250°C. The volatile thermal degradation products yielded by the chemical breakdown of carbon fibre based epoxy composite at the glass transition temperature were measured, where a mixture of aromatic compounds and unsaturated compounds were identified and quantified.

Conclusions

The sensors have demonstrated a good performance to real-time detection of gases, although cross-response of sensors to a range of volatiles might occur. It is important to accurately determine the characteristics of these compounds in aircraft cabins, and modelling is needed to link the emission rate in such tests with a possible concentration in the cabin.

Novel Aspect

A methodology combining real-time data information and GC-MS data was set to investigate further on potential emissions of volatiles. To the authors knowledge, this work has not been previously reported.

References

Dastin S.J. Aircraft Composite Materials and Structures. In: Fridlyander I.N., Kostikov V.I. (eds) MICC 90. Springer, Dordrecht (1991).

MOr.06 - Clinical chemistry

Chairs: Hugo Rocha, Giancarlo la Marca

MOr.06-1 **Keynote: TANDEM MASS SPECTROMETRY IN THE CLINICAL LABORATORY: CAN STANDARDIZATION IMPROVE OPERATION**

Marzia Pasquali

University of Utah School of Medicine

Keywords: Mass Spectrometry, Standardization, Operational Efficiency

Today's clinical laboratory is focused on the delivery of accurate results at the right time to the right user. To reach this goal, all components of the laboratory, from the technologists on the bench to the analyst handling the LIS (laboratory information system) need to work in an integrated and coordinate way. The large number of methodologies used in the clinical

laboratory can make the integration challenging and standardization processes are becoming the answer to the challenges brought by this diversity. Among the methodologies, mass spectrometry, whether single stage or tandem, often coupled with liquid or gas chromatography, is becoming one of the most widespread technique in the clinical laboratory, because of its sensitivity, specificity, rapid turnaround time, and low sample volume requirement. Mass spectrometry is therefore used in many areas, from newborn screening to therapeutic drug monitoring, endocrinology, infectious diseases, in addition to being used as a research tool in metabolomics applications.

We have evaluated standardization processes for the mass spectrometry area in our laboratory, in which we have two opposite settings: 1) newborn screening, high throughput/rapid turnaround time/low sample volume and 2) biochemical genetics, low throughput/high complexity/rapid turnaround time.

From our preliminary observations, we have seen that standardization of the platform used may be beneficial in the first setting, while in the second setting, standardized analytical/instrument methods could improve the workflow, independently from the platform.

In both cases, automated data transfer to the LIS is a key factor in delivering test results rapidly and accurately.

The implementation of standardization should take into account the laboratory setting in order to optimize its efficiency. Each laboratory should identify “why standardize?” and “what to standardize?” before identifying “how to standardize”. This will fulfill the mission of the clinical laboratory by improving operational efficiency without compromising delivery of care.

MOr.06-2 A COMBINED MICRO-EXTRACTION BY PACKED SORBENT-UHPLC-ESI-MS/MS PROCEDURE FOR THE ANALYSIS OF PROSTANOIDS AND ISOPROSTANOIDS IN DRIED BLOOD SPOTS OF PRETERM NEWBORNS WITH PATENT DUCTUS ARTERIOSUS

Denise Biagini⁽¹⁾ - ***Shaula Antoni***⁽¹⁾ - ***Silvia Ghimenti***⁽¹⁾ - ***Tommaso Lomonaco***⁽¹⁾ - ***Francesca Giuseppa Bellagambi***⁽¹⁾ - ***Armando Cuttano***⁽²⁾ - ***Rosa Teresa Scaramuzzo***⁽²⁾ - ***Massimiliano Ciantelli***⁽²⁾ - ***Fabio Di Francesco***⁽¹⁾ - ***Roger Fuoco***⁽¹⁾

⁽¹⁾ University of Pisa, Department of Chemistry and Industrial Chemistry, Pisa - ⁽²⁾ AOUP Santa Chiara Hospital, University of Pisa, Division of Neonatology, Pisa

Keywords: tandem mass spectrometry, micro-extraction by packed sorbent, dried blood spots, isoprostanoids, preterm newborns

Introduction

An innovative UHPLC-ESI-MS/MS method for the analysis of prostaglandin-E₂, 8-iso-prostaglandin-F_{2α} and 8-iso-prostaglandin-E₂ in dried blood spots is proposed. The non-invasive and painless monitoring of these inflammation and oxidative stress mediators in preterm newborns suffering from Patent Ductus Arteriosus could assess their potential role as predictors of response to therapy [1].

Methods

Analytes were extracted from dried blood spots (20 µL x 3) with 70:30 v/v methanol:water mixture. The extracts were purified by a semi-automated micro-extraction by packed sorbent (MEPS) before injection in the UHPLC-MS/MS instrument. Analytes were separated by a reversed-phase Polaris-C18 column and a gradient elution with aqueous formic acid (0.1%) and 50:50 v/v methanol:acetonitrile mixture. Triple quadrupole mass spectrometer operated in ESI(-) and MRM mode.

Results

The fast and easy sample clean-up by MEPS guaranteed a high reproducibility (RSD < 5%) and extraction yield (>70%) in less than ten minutes, and limited the required volume of both

solvents (30-100 μL) and sample (100 μL), with respect to conventional solid phase extraction [2]. The MRM acquisition ensured high sensitivity and selectivity and the presence of the analytes in the samples was confirmed by a q/Q ratio within $\pm 20\%$ of the value measured in reference standards. Limits of detection below 20 pg/mL , linearity ($R^2 > 0.99$) over three orders of magnitude and satisfactory overall intra- and inter-day precisions (RSD $< 10\%$) were obtained. The deuterated internal standard was successfully laid on the filter paper instead of the common addition to the extraction solvent [2], in order to control both storage and extraction steps [3]. Finally, particular attention was paid to stability issues.

Conclusions

To the best of our knowledge, this is the first analytical method combining a fast and reproducible MEPS procedure and a highly sensitive and selective UHPLC-ESI-MS/MS analysis for the determination of prostaglandin- E_2 , 8-iso-prostaglandin- $\text{F}_{2\alpha}$ and 8-iso-prostaglandin- E_2 in dried blood spots. The method is suitable for monitoring preterm newborns with Patent Ductus Arteriosus.

Novel Aspect

The very reliable combination of tandem mass spectrometry and fast, semi-automated, reusable MEPS for prostanoids and isoprostanoids analysis in dried blood spots.

References

1. Inayat M., Bany-Mohammed F., Valencia A., Tay C., Jacinto J., Aranda J.V., Beharry K.D., *Am J Perinatol.*, 32, 895-904 (2015).
2. Bastani N.E., Gundersen T.E., Blomhoff R., *Rapid Commun. Mass Spectrom.*, 26, 645-652 (2012).
3. Abu-Rabie P., Denniff P., Spooner N., Brynjolfsson J., Galluzzo P., Sanders G., *Anal. Chem.*, 83, 8779-8786 (2011).

MO.06-3 PROFILING OF MEMBRANE LIPIDS RELEVANT TO ALZHEIMER'S DISEASE PATHOGENESIS

Zdenek Spacil⁽¹⁾ - **Eliska Cechova**⁽¹⁾ - **Petr Telensky**⁽²⁾ - **Seweryn Olkowicz**⁽²⁾ - **Lukas Opalka**⁽³⁾ - **Katerina Sheardova**⁽⁴⁾ - **Radim Jancalek**⁽⁵⁾ - **Jana Klánová**⁽¹⁾ - **Ales Hampl**⁽⁵⁾ - **Jiri Damborsky**⁽¹⁾

⁽¹⁾ Masaryk University, Faculty of Science, Research Centre for Toxic Compounds in the Environment, Brno – ⁽²⁾ Charles University, Faculty of Science, Prague - ⁽³⁾ Charles University, Faculty of Pharmacy, Inorganic and Organic Chemistry, Hradec Kralove - ⁽⁴⁾ Saint Ann's Faculty Hospital, International Centre for Clinical Research, Brno - ⁽⁵⁾ Masaryk University, Faculty of Medicine, Brno

Keywords: metabolic profiling, gangliosides, Alzheimer's disease, monosialoganglioside GM1, membrane lipidomics

Introduction

Alzheimer's disease (AD) is the most common cause of dementia (60-80% of cases) in elderly population worldwide [1]. The current body of research supports the concept that protein aggregation initiates the onset of AD [2]. However, in spite of the large number of studies, a mechanism of pathogenesis and potential treatment for AD remains elusive. We explored the distribution of membrane lipids (gangliosides) as they supposedly trigger protein aggregation [3].

Methods

We applied ultra-high performance liquid chromatography (UHPLC) and selected reaction monitoring (SRM) as well as high resolution/accurate mass (HR/AM) orbital ion trap mass spectrometry to characterize membrane lipid composition (GM1, GM2, GM3, GD1a/b, GD2, GD3, GT1b, GQ1b) in human clinical samples (i.e. serum, cerebrospinal fluid - CSF, brain tissue)

and samples from animal model of aging (*Sprague Dawley, Wistar*). We compared 1 month and 11 months old animals (n=26).

Results

Our results show optimization of ganglioside extraction protocol and development of UHPLC-SRM assay. We tested several solvent systems to achieve approx. 70% extraction efficacy in biological matrix. UHPLC was optimized and suitable analytical column and mobile phase additives chosen to allow for separation of mono-, di-, tri- and quaterosialo species. Gangliosides present in biological samples were identified using HR/AM mass spectrometry. Consequently positive and negative ion mode SRM assays were utilized for determination of gangliosides. The dominant brain gangliosides were GM1, GD1a/b and GT1b and similar composition was found in CSF. In serum the most abundant species are GM3, GM2, GD1a/b, GD2 and GT1b. We have used custom synthesized isotopically labelled internal standards for absolute quantification in biological samples.

Conclusions

New analytical method established to determine 9 major gangliosides in brain tissue, CSF and serum samples with nanomolar LODs. The differences in ganglioside composition were demonstrated in specific brain structures (i.e. cerebral cortex, cerebellum, hypothalamus, hippocampus, brain stem). The composition of gangliosides in brains reflected CSF and partially also serum concentrations. The age dependent changes in composition of major gangliosides were evaluated.

Novel Aspect

Novel analytical approach to study distribution of important membrane lipids, which are supposedly playing role in Alzheimer's disease pathogenesis.

References

1. Mayeux, R. & Stern, Y. *Cold Spring Harb. Perspect. Med.* **2**, a006239–a006239 (2012).
2. Koo, E. H., Lansbury, P. T. & Kelly, J. W. *Proc. Natl. Acad. Sci.* **96**, 9989–9990 (1999).
3. Cebecauer, M., Hof, M. & Amaro, M. *Biophys. J.* **113**, 1194–1199 (2017).

MOr.06-4 THE USE OF STABLE ISOTOPES IN THE STUDY OF ALBUMIN METABOLISM IN PATIENTS WITH ACUTE LUNG DISEASE

Manuela Simonato⁽¹⁾ - **Sonia Giambelluca**⁽²⁾ - **Luca Vedovelli**⁽¹⁾ - **Lorenza Dalla Massara**⁽³⁾ - **Carlo Ori**⁽³⁾ - **Virgilio Paolo Carnielli**⁽⁴⁾ - **Paola Cogo**⁽⁵⁾

⁽¹⁾ Fondazione Istituto di Ricerca Pediatrica Città della Speranza, PCare Laboratory, Padova - ⁽²⁾ Padova University Hospital, Department of Women's and Children's Health, Padova - ⁽³⁾ Padova University Hospital, Department of Medicine, Section of Anesthesiology and Intensive Care, Padova - ⁽⁴⁾ Polytechnic University of Marche and Azienda Ospedaliero-Universitaria Ospedali Riuniti, Division of Neonatology, Department of Clinical Sciences, Ancona - ⁽⁵⁾ University of Udine, Department of Medicine, Udine

Keywords: acute respiratory distress syndrome; albumin; isotope ratio mass spectrometry; stable isotope

Introduction

Acute respiratory distress syndrome (ARDS) is a severe form of lung injury characterized by an increased inflammatory burden and changes in the alveolar-capillary membrane permeability with extravasation of fluid and protein into the alveoli[1]. Despite the clear association between plasmatic hypoalbuminemia and poor outcome, there are no data on albumin metabolism in ARDS patients.

Methods

^{13}C -Leucine, as albumin precursor, was administered to 13 ARDS adults and to 8 patients without lung disease (controls). Albumin was isolated from serial plasma and tracheal aspirate (TA) samples, and its kinetic parameters were derived from the ^{13}C enrichment curve obtained by isotope ratio mass spectrometry[2]. As index of alveolar-capillary permeability the TA/plasma albumin ratio was measured.

Results

In patients with ARDS the fractional synthetic rate was significantly higher than in control patients [15.1 (12.3-22.7) vs 11.9 (8.9-12.9) %/day; $p=0.04$].

Median albumin secretion time and half-life were 0.7 (0.3-1.9) vs 0.6 (0.3-1.0) h, and 43.3 (35.7-64.9) vs 42.7 (33.3-76.0) h in ARDS and in control patients respectively, no significant differences were observed.

In patients with ARDS the TA to plasma albumin ratio measured in the early phase after the diagnosis of acute lung injury was 7-fold times higher the level found in control patients ($p<0.001$). In the ARDS group the TA/plasma albumin ratio was significantly correlated with the $\text{PaO}_2/\text{FiO}_2$, as index of lung dysfunction ($r=-0.732$, $p=0.016$).

Conclusions

In this study, stable isotopes have been used to better understand albumin metabolism during a disease state like ARDS. Albumin synthesis rate in ARDS patients was higher than in control adults. Furthermore, the comparison between TA/plasma albumin ratio of ARDS and control patients revealed a dramatic increase in alveolar-capillary membrane permeability during the genesis of ARDS.

Novel Aspect

This study will give new insight into albumin turnover and compartmental distribution in patients with ARDS.

References

1. Ware L.B., Matthay M.A., *The New England Journal of Medicine*, 342, 1334 (2000).
2. Simonato M., Baritussio A., Pioselli B., Ori C., Catinella S., Carnielli V.P., Cogo P., *Analytical and Bioanalytical Chemistry*, 406, 6225 (2014).

MOr.06-5 MASS-SPECTROMETRY APPROACH FOR RAPID INTRAOPERATIVE IDENTIFICATION OF BRAIN TUMOR TISSUES - A WAY FOR ONLINE ANALYSIS

Stanislav Pekov⁽¹⁾ - **Vasily Eliferov**⁽¹⁾ - **Evgeny Zhvansky**⁽¹⁾ - **Alexander Vorobyev**⁽¹⁾ - **Vsevolod Shurkhay**⁽²⁾ - **Anatoly Sorokin**⁽¹⁾ - **Igor Popov**⁽¹⁾ - **Alexander Potapov**⁽²⁾ - **Evgeny Nikolaev**⁽³⁾

⁽¹⁾ Moscow Institute of Physics and Technology, Dolgoprudnyj - ⁽²⁾ N. N. Burdenko Scientific Research Neurosurgery Institute, Moscow - ⁽³⁾ Skolkovo institute of science and technology, Skolkovo

Keywords: ambient mass-spectrometry, clinical application, brain tumor, neurosurgical support, microextraction.

Introduction

It is known that the lipid metabolism reprogramming is one of the hallmarks of cancer so mass-spectrometry-based method of lipid profiling is a promising approach for precise intraoperative tumor identification. In this research we propose the online method of tissue profiling based on

express and simple capsular extraction followed by ESI with an identification algorithms based on a search through a reference database containing mass spectra of different tumors.

Methods

The online microextraction system consists of disposable capsules compatible with any standard ESI source. Tissues, dissected by the surgeon were placed into the capsules which are then connected to the ESI interface for analysis. Features were extracted from mass-spectrometry data and aligned for further analysis. Significant features were detected by shrinkage discriminant analysis and classifiers were developed with lasso-regression.

Results

It was demonstrated that it is possible to prepare large amount (up to 5 per minute) of capsules with samples immediately during the operation because of simplicity of placing sample into the preassembled capsules so almost all the delay between receiving the tissue from surgeon and analysis depends only from transportation time between operational suite and laboratory room. The minimum amount of sample is very small (less than 1 cubic mm) so proposed method could be used during neurosurgery intervention where amount of resection could be very small. The data obtained with proposed microextraction method are in good correlation with our previous results obtained with spray-from-tissue ionization source. The rapid identification of the tumor type is allowed by the search inside the collected reference database containing mass spectra of samples that are also characterized by histopathological methods.

Conclusions

Application of microextraction capsules significantly simplifies the manipulations with samples by excluding any special sample preparation steps. Rapidity of method is also very important during the surgery as the minimization of surgery duration improves patient's postoperative prognosis. Developed data processing algorithms and classifiers provide the reliable differentiation of extracted brain tissue by their lipid profiles.

Novel Aspect

The online mass-spectrometry profiling of brain tissues could be used as the intraoperative technology for characterization and identification of removed tissues.

MOr.07 - Quantitative proteomics

Chairs: Michael Linscheid, Nathalie Norais

MOr.07-1 **Keynote: DRUG EFFECTS ON PROTEIN HOMEOSTASIS**

Marcus Bantscheff

GlaxoSmithKline, GSK Cellzome GmbH, Heidelberg

Keywords: protein turnover, drug effects, protein homeostasis, targeted degraders

Introduction

Protein degradation plays important roles in biological processes and is tightly regulated. Further, targeted proteolysis is an emerging research tool and therapeutic strategy. However, proteome-wide technologies to investigate the causes and consequences of protein degradation in biological systems are lacking.

Methods

mPDP enables the sensitive and comprehensive relative quantification of cellular mature and nascent protein pools across a range of different treatment conditions in biological duplicates and in a single mass spectrometric experiment by combining dynamic SILAC with chemical labeling using neutron-encoded tandem mass tags (TMT).

Results

We developed ‘multiplexed proteome dynamics profiling’, mPDP [1], a mass spectrometry-based approach combining dynamic-SILAC labelling with isobaric mass tagging for multiplexed analysis of protein degradation and synthesis. When applied combination with other quantitative proteomics approaches such as chemoproteomics and thermal protein profiling, mPDP provides unique insights in drug mechanism-of-action. The method was applied in three proof-of-concept studies, and uncovered different responses induced by the bromodomain inhibitor JQ1 versus a JQ1-proteolysis targeting chimera. Further, we elucidated distinct modes-of-action of estrogen receptor modulators, and we comprehensively classified HSP90 clients based on their requirement for HSP90 constitutively or during synthesis demonstrating that constitutive HSP90 clients have lower thermal stability than non-clients, higher affinity for the chaperone, vary between cell types and change upon external stimuli.

Conclusions

We have established multiplexed proteome dynamics profiling a powerful means to elucidate the mechanism of action of bioactive molecules, thus complementing current proteome-wide target identification strategies. mPDP identifies targets of protein degraders and more generally enables the discovery of regulatory degradation mechanisms in biological systems.

Novel Aspect

PROTAC target identification and detection of regulated degradation events.

References

Savitski MM et al., *Cell* **173**, 260-274 (2018)

MOr.07-2 MEASURING NASCENT PROTEOMICS IN SIGNALING AND CELLULAR DIFFERENTIATION

Alma Burlingame⁽¹⁾ - **Nancy Phillips**⁽²⁾ - **Craig Forester**⁽³⁾

⁽¹⁾ UCSF, Pharmaceutical Chemistry, San Francisco - ⁽²⁾ UCSF, Pharmaceutical Chemistry, San Francisco - ⁽³⁾ UCSF, Urology, San Francisco

Keywords: *O*-propargylpuromycin, nascent proteome, erythropoiesis, mass spectrometry

Introduction

Regulation of gene expression is controlled by mRNA translation and thus defines the functional proteome. Experimental methodologies are needed that can rapidly capture and reveal the actual nascent proteomes. In addition, such methods must be able to measure changes in nascent protein abundances quantitatively in response to pharmacologic or environmental cues.

Methods

We have developed a protocol using an analog of puromycin, a tyrosine-tRNA mimetic, called *O*-propargylpuromycin (OPP), to identify nascent proteomes rapidly. OPP is cell permeable and catalytically incorporated into nascent polypeptide chains as a specific covalent label during elongation resulting in premature termination. These polypeptides are biotinylated using click chemistry, then isolated and measured by tandem mass spectrometry[1].

Results

We have chosen the erythropoiesis paradigm known to be dependent upon an intricate balance involving maintenance of erythroid precursors and their rapid differentiation to accommodate hematopoietic stresses. In this study, we employed extensively self-renewing erythroblast (ESRE) cells to define the earliest proteomic program unique to the decision between expanding erythroid precursors and terminal maturation. Using OPP-ID, we discovered a network of genes involved in splicing, nucleic acid secondary structure unwinding, and cell cycle regulation which are rapidly repressed in the absence of dexamethasone. While we observed genes previously known to be associated with early hematopoietic precursor states (such as Ddx18, Stag2 and Pold1), we also identified multiple genes such as Prpf6, Sart1, Psen2, and Mycbp2 which have not been previously implicated in erythropoiesis.

Conclusions

This OPP-ID strategy provides rapid experimental view of changes in gene expression at the nascent protein level. This opens a new window for the determination of the actual proteomic changes occurring in response to a variety of cues that direct cellular fates such as proliferation or differentiation.

Novel aspect

This OPP-ID strategy has the potential to ask innumerable questions on the composition and rapid dynamics of proteomic networks at the organismal level.

References

1. Forester c., et al., *Proc. Natl. Acad. Sci.*, 2018, doi: 10.1073/pnas.1707514115

MOr.07-3 **A CHEMICAL PROBE FOR UNDERSTANDING HOW PROTEOME FOLDEDNESS CHANGES UNDER PROTEOSTASIS STRESS**

Gavin Reid⁽¹⁾ - Dezerae Cox⁽²⁾ - Yuning Hong⁽³⁾ - Danny Hatters⁽²⁾

⁽¹⁾ University of Melbourne, School of Chemistry, Department of Biochemistry and Molecular Biology, Parkville - ⁽²⁾ University of Melbourne, Department of Biochemistry and Molecular Biology, Parkville - ⁽³⁾ La Trobe University, Department of Chemistry and Physics, Bundoora

Keywords: proteostasis, protein folding, chemical proteomics

Introduction

The components of the cellular quality control network responsible for maintaining proteostasis, some of which function to protect the folded state of proteins, are now well understood. However, it is still unknown which proteins in the proteome are most vulnerable to unfolding, misfolding, or aggregation, in response to diverse proteostasis stresses.

Methods

Here, we demonstrate the application of a thiol-specific probe, tetraphenylethene maleimide (TPE-MI), coupled with stable isotope labelling and MS-based quantitative proteomic analysis, for measuring changes in the reactivity of free cysteine residues that are normally buried in the core of globular proteins but become exposed upon unfolding due to proteostasis stress¹.

Results

We have demonstrated the potential of this method to (i) probe the kinetics of single-protein folding/unfolding for endogenous proteins on a proteome-wide scale, and (ii) monitor the specific changes in cellular proteome foldedness that occur in response to a variety of stressors (e.g., ER stress (tunicamycin), Heat Shock, HSP-70 inhibition, proteasome inhibition, or

overexpression of a non-native misfolded protein, barnase²). Notably, we show that different stress conditions result in little overall similarity in changes in foldedness, and share only as small subset of common vulnerable core machinery.

Conclusions

TPE-MI is demonstrated as a promising tool to probe proteostasis at both the proteome wide and individual protein scales. This knowledge will contribute to our understanding of disorders characterised by proteostasis imbalance, and will assist in targeting those proteins most prone to misfolding under stress conditions.

Novel Aspect

This approach enables analysis of folding/unfolding kinetics for endogenous proteins on a proteome-wide scale, without expression of specific protein reporters.

References

1. Chen, M., Moily, N. S., Bridgford, J., Wood, R., Radwan, M., Smith, T. A., Song, Z., Tang, B. Z., Tilley, L., Xu, X., Reid, G. E., Pouladi, M. A., Hong, Y., Hatters, D. M. *Nature Communications*, 8, 474 (2017).
2. Wood, R.J., Ormsby, A.R., Radwan, M., Cox, D., Sharma, A., Vöpel, T., Ebbinghaus, S., Oliveberg, M., Reid, G.E., Dickson, A., Hatters, D.M. *Nature Communications*, 9, 287 (2018).

MOr.07-4 **TMT QUANTIFICATION IN MS2 USING HIGH FIELD ASYMMETRIC WAVEFORM ION MOBILITY SPECTROMETRY**

Sibylle Pfammatter - Bonneil Eric - Pierre Thibault

Institute for Research in Immunology and Cancer (IRIC), Université de Montréal, Montréal

Keywords: Proteomics, FAIMS, PTM, TMT quantification

Introduction

Proteomic analyses are often limited by the overwhelming proportion of confounding background ions that compromise the identification and quantification of low abundance peptides. Gas phase ion fractionation using high field asymmetric waveform ion mobility (FAIMS) can reduce sample complexity [1], and improve quantitative proteomic measurements [2]. Here we report a new FAIMS interface that further expand the sensitivity and depth of proteomic analyses.

Methods

LC-MS/MS experiments were performed with and without FAIMS on an Orbitrap Fusion mass spectrometer. We used cylindrical FAIMS electrodes (1.5 mm gap) with both electrodes set to 100°C, pure nitrogen as carrier gas and a dispersion voltage (DV) of -5000 V. LC-FAIMS-MS/MS experiments were conducted by stepping through 2-4 different CVs in a single run. Cell extracts were reduced, alkylated and digested with trypsin. TMT at 1:2 label excess was done in HEPES buffer, pH 8.3.

Results

To evaluate the advantages of FAIMS, we first compared the number of identification obtained for triplicate injections of a 500 ng HEK293 digest analyzed by LC-MS/MS with and without FAIMS. These analyses revealed that FAIMS provided a 30% enhancement in identification with 39791 peptides (2876 protein groups) compared to 31397 peptides (2414 protein groups) without FAIMS. We next used isobaric peptide labeling with TMT to profile the dynamic changes in protein abundance for HEK293 cells exposed to heat shock. A different TMT label was used for

each 1h interval for a treatment period up to 9h. Samples were analyzed using both synchronous precursor selection (SPS-MS3) [3] and LC-FAIMS-MS/MS. SPS-MS3 enabled the quantification of 12400 peptides corresponding to 1229 proteins, of which 375 proteins showed dynamic changes in abundance. FAIMS further enhanced the comprehensiveness of quantitative proteomics by more than 2-fold with 30848 peptides assigned to 2646 proteins, of which 902 proteins showed dynamic changes in abundance.

Conclusions

A new FAIMS interface was coupled to the latest generation of Orbitrap instruments and enabled unparalleled sensitivity in proteomics analyses where multiple CVs can be combined in a single LC-MS/MS run. The enhancement in ion transmission and duty cycle of the extended the depth of proteomic analyses by at least 30% and expand the comprehensiveness of multiplex quantitative proteomics using isobaric peptide labeling by at least 2-fold compared to SPS-MS3 analyses.

Novel Aspect

Novel FAIMS interface combined with an Orbitrap Tribrid MS improves the sensitivity of proteome analyses and the comprehensiveness of quantitative proteomic measurements.

References

1. Bonneil E., Pfammatter S., Thibault P., *Journal of Mass Spectrometry*, 15, 4653-65 (2015)
2. Pfammatter S., Bonneil E., Thibault P., *Journal of Proteome Research*, 50, 1181-95 (2016)
3. McAlister et al., *Anal. Chem.*, 86, 7150–7158 (2014)

MOr.07-5 THE KEY TO BIG BIOBANKS: ANALYZING 24 PROTEOMES PER DAY BY MICRO-FLOW SWATH® ACQUISITION AND SPECTRONAUT PULSAR ANALYSIS

Jan Muntel⁽¹⁾ - **Nick Morrice**⁽²⁾ - **Roland Bruderer**⁽¹⁾ - **Lukas Reiter**⁽¹⁾

⁽¹⁾ Biognosys AG, R&D, Schlieren - ⁽²⁾ AB SCIEX UK Limited, AB SCIEX, Warrington

Keywords: DIA, SWATH, high throughput proteomics, micro-flow, FFPE

Introduction

Micro-flow chromatography provides a similar workflow sensitivity compared to nano-flow separations but greatly increased robustness enabling high throughput proteomics for 100's of samples. In combination with data-independent acquisition methods (DIA) such as SWATH® acquisition, 1,000s of proteins can be quantified with high data completeness and high quantitative precision. Here, we show the potential of the workflow for 105 cancer samples.

Methods

FFPE colon tissue samples (10 healthy, 95 cancer) were obtained from a public repository. Proteins were extracted as described previously [1]. Tryptic digests were spiked with iRT peptides (Biognosys) and separated on a NanoLC 425 system (SCIEX) by a 43 min gradient at a flow rate of 5µl/min coupled to a TripleTOF® 6600 system (SCIEX) operated in SWATH mode. Run time per sample was 1h. Data were analyzed in Spectronaut Pulsar (Biognosys) with a project specific library.

Results

The sample preparation resulted in average in >100 µg peptides per 10 µm slice, which is plenty for the micro-flow SWATH® acquisition workflow. The whole dataset (105 samples) was acquired in only 4.5 days. Data analysis in Spectronaut Pulsar took one additional day demonstrating the

large potential of the workflow for even more samples. In total, about 4,500 protein groups were quantified at 1% peptide and protein FDR (3,600 in average per run) providing a sufficient depth for biological interpretation of the data. Data analysis revealed that a larger number of proteins were expressed in the cancer cohort. Unsupervised clustering of the data separated the healthy and the cancer cohort. Interestingly two samples of the healthy cohort clustered together with the cancer cohort. Within the cancer cohort, three major groups were identified, potentially different cancer subtypes. The data will be further biologically interpreted using ingenuity pathway analysis (IPA).

Conclusions

This study demonstrates, how high throughput proteomics can be used to analyze large sample sets from tissue biobanks. Application of the micro-flow SWATH® acquisition workflow combined with Spectronaut Pulsar data analysis generates data from these large resources in a timely fashion and therefore pave the way to a better understanding of health and disease.

Novel Aspect

Micro-flow/SWATH® acquisition and Spectronaut Pulsar analysis enables efficient and comprehensive analysis of large sample cohorts (>100 samples).

References

1. Buczak K., et al., *Molecular & Cellular Proteomics* (2018)

MOr.08 - Chiral discrimination

Chairs: Anne Zehnacker-Rentien, Maurizio Speranza

MOr.08-1 **Keynote: CHIRAL RECOGNITION IN PROTONATED MOLECULAR CLUSTERS STUDIED BY MASS-COUPLED IR-UV DOUBLE RESONANCE SPECTROSCOPY WITH COLD ION TRAP**

Masaaki Fujii

Laboratory for Chemistry and Life Science, Institute of Innovative Research, Tokyo Institute of Technology

Keywords: Chirality, Molecular Recognition, ion trap, IR, laser spectroscopy

Introduction

Understanding stereochemical effects on the structure of molecular aggregates provides insight into the role of chirality in shaping biological supramolecular systems at the molecular level. Here we focus on 1-amino-2-indanol (AI2H^+), a proto- typical cyclic 1,2-amino-alcohol, and reveal how the the chirality of the monomer subunits controls the overall structure of its dimers [1].

Methods

Here, we apply IR-UV double-resonance laser spectroscopy coupled with ESI / cold ion trap mass spectrometry to the protonated homodimers of (1R,2S)-cis- and (1R,2R)-trans- AI2H^+ (c- AI2H^+ , t- AI2H^+). The IR-UV double-resonance spectroscopy enables us to measure UV and IR spectra of a single conformer. Such conformer selection is a big advantage of mass-coupled laser spectroscopy [1].

Results

The conformer-selected UV spectra are measured for both cis- and trans-monomers (c- AIH^+ , t- AIH^+) by IR-UV hole burning spectroscopy. The structural analysis with help of quantum chemical

calculations at B3LYP-D3/6-31G++(d,p) level show that c-AIH⁺ has an intramolecular NH⁺...O hydrogen bond (H-bond), while t-AIH⁺ lacks such an interaction. Conformer-selected UV spectra of cis- and trans-AIH⁺ homodimers are measured and three and two conformers are assigned, respectively. Their structures are analyzed by the comparison between conformer selected IR spectra and theoretical IR spectra. The obtained structures of cis- and trans-dimers are clearly different. The structural difference is rationalized by the presence and absence of the intramolecular H-bond in cis and trans monomer subunits.

Conclusions

This work demonstrates that the combination of ESI cold ion trap mass spectrometry and laser spectroscopy is a powerful tool to study the structure of molecular aggregates and effect of chirality. The methods in this work will allow us to expand the size of molecular systems to biological supramolecular systems. In the presentation, the application to functional peptides and its complex with chiral ligands will also be shown [2].

Novel Aspect

The approach developed in this work will provide us a deep understanding of the stereochemistry-based control of biological supramolecular system.

References

1. Bouchet A. et al., *Phys. Chem. Chem. Phys.*, published online, DOI: 10.1039/c8cp00787j.
2. Sekiguchi T. et al., *Ang. Chem. Int. Ed.* published online, DOI: 10.1002/anie.201712976.

MOr.08-2 A NEW MS-BASED HYPHENATED METHOD FOR ANALYSIS OF HEXURONIC ACID EPIMERS AND SULFATE PATTERNS IN GLYCOSAMINOGLYCAN OLIGOSACCHARIDES BY MSⁿ COUPLED TO INFRARED ION SPECTROSCOPY

Gina Renois-Predelus⁽¹⁾ – ***Baptiste Schindler***⁽¹⁾ – ***Loïc Barnes***⁽¹⁾ – ***Stéphane Chambert***⁽²⁾ – ***Abdul-Rahman Allouche***⁽¹⁾ – ***Isabelle Compagnon***⁽¹⁾

⁽¹⁾ *Institut Lumière Matière, Université Claude Bernard Lyon 1, Lyon* - ⁽²⁾ *Institut de Chimie et de Biochimie Moléculaires et Supramoléculaires, Université Claude Bernard Lyon 1, Lyon*

Keywords: Glycosaminoglycan, MS/IR, hexuronic acid epimers, structural characterization, sulfate patterns

Introduction

We report an original MS-based hyphenated method for the elucidation of sulfate position and the epimerization of hexuronic acids in glycosaminoglycan (GAG) fragments.

Methods Coupling of mass spectrometry with IRMPD spectroscopy

This method consists of measuring simultaneously the MS/MS spectrum and the gas phase IR spectrum to gain direct structural information. This is possible using a customized MS instrument, modified to allow injection of a tunable IR laser inside the instrument for in situ spectroscopy of trapped ions.

Results

This approach was tested for the analysis of the hexuronic fragments composing the hyaluronic acid tetrasaccharide [1] and the sulfated GalNAc fragments present in the chondroitin sulfate disaccharides and dermatan sulfate tetrasaccharide [2].

Conclusions

We show that the gas phase IR fingerprint of reference hexuronic acid and sulfated GalNAc monosaccharides are distinct and can be used to detect them in oligosaccharides.

Novel Aspect

The method presented here is predictive and allows structural elucidation of unknown GAG fragments, even in the absence of referenced standards.

References

1. Schindler B., Renois-Predelus G., Bagdadi N., Melizi S., Barnes L., Chambert S., Allouche A.-R., Compagnon I., *Glycoconjugate Journal*, 34, 421–425 (2017).
2. Renois-Predelus G., Schindler B., Compagnon I., *J. Am. Soc. Mass Spectrom.* (Publication in press 2018)

MO.r.08-3 **CHIRAL DISCRIMINATION IN COLLISION-INDUCED DISSOCIATION EXPERIMENTS: A THEORETICAL APPROACH OF THE CYCLOPHE-HISH⁺ SYSTEM**

Ariel Pérez-Mellor - Riccardo Spezia

Laboratoire de Chimie Théorique, LCT, Sorbonne Université, CNRS, F. 75005 Paris, France

Keywords: Collision-induced dissociation (CID), diketopiperazine (DKP), chiral discrimination, Born-Oppenheimer molecular dynamic (BOMD).

Introduction

This theoretical study is inspired by the strong experimental differences observed in the breakdown curves of the collision-induced dissociation of the two diastereoisomers of the diketopiperazine peptide cycloSPhe-SHisH⁺ (SS) and cycloSPhe-RHisH⁺ (SR) (unpublished results). The major chiral discrimination is observed in the carbon monoxide and ammonia loss channels. For the former case, the fragmentation efficiency of the SR diastereomer is larger in contrast to what happens in the latter.

Methods

A comprehensive picture of unimolecular dissociation is obtained by a Born-Oppenheimer molecular dynamic simulation. Internal energy activation through a micro-canonical ensemble and the RM1-D semi-empirical Hamiltonian has been used for this purpose. The method used showed good a balance between accuracy and computational time. All the simulations were carried out coupling VENUS (for dynamics) and MOPAC (for electronic structure calculation) packages.

Results

Within this level of theory, we are able to reproduce all the main fragmentation channels observed in the CID experiment starting from the geometry determined by infrared multiple photon dissociation spectroscopy. This geometry agrees with the calculated global minimum for both systems and shows that the proton lies in the histidine residue. The main fragmentation products are loss of ammonia, loss of water, loss of CO, loss of formamide or a combination of NH₃+CO (these two last products have the same mass), loss of carbamic acid or a combination of isocyanic acid + H₂O, loss of HCONH₂+CO or a combination of NH₃+2CO, the formation of the iminium ion of Phe and the formation of the iminium ion of His. The Histidine residue plays an important role. It consists in carrying the proton to highly energetic proton sites like oxygen and nitrogen atoms in DKP ring (following the proton mobile model). There is a strong difference in the mobility of the proton in the two diastereomers due to the chirality constrain. In the SS case, the proton mobility is lower in comparison to SR case due to the fact the two residues remain on the same side of the DKP ring interacting. Then, the SR system shows, overall, higher fragmentation yield than SS system in agreement with experiments. This is likely due to a lower activation energy.

Conclusions

The strong experimental differences observed in the CID experiments can be rationalized thanks to theoretical calculations which point out the role of the Histidine residue in the proton mobility which is related to differences between the two diastereoisomers.

Novel Aspect

This chiral discrimination has been noticed before but a clear rationalization was not found. Calculations can help in understanding it from a molecular basis.

MOr.08-4 **ENANTIOSELECTIVE REDUCTION FROM QUATERNARY COPPER D,L AMINO ACID COMPLEXES: ORIGIN OF THE WEAK CHIRAL EFFECT OBSERVED WITH D/L ASPARTIC AND GLUTAMIC ACIDS**

Jean-Claude Tabet

Keywords: chirality, enantioselective reduction, CID, interpretation

Introduction

Amino acid enantiomers (AA_L/AA_D) distinction was a challenge solve after the 2000s. A method proposed by G. Cooks (1) was based on (i) electrospray formation of quaternary $[2AA'_L+AA_{L/D}-H+Cu^{(II)}]^+$ complex diastereomers (AA'_L as reference, $AA_{L/D}$, that of unknown chirality) and (ii) their comparative dissociations by CID. Applying this approach and using AA'_L =phenylglycine, an enantioselective reduction accompanying expected chiral effects were observed (2).

Methods

Each AA was dissolved with $CuCl_2$ in MeOH/ H_2O (50:50) within 20 pmol.mL⁻¹. HD exchanges were performed. Experiments were done positive ESI with ion trap, and with LTQ/Orbitrap. Solutions were respectively infused at a flow rate of 160 mL/h and 100 μ L/h. Resonant excitation (CID) and non-resonant (HCD) were used. The source settings are as follows: capillary voltage, -3.5 kV; heat temperature, 100°C; The image signal was amplified and digitized using 1M data points.

Results

Using the Cooks' approach (1), the phenylglycine (Phg_L)(2) was used as chiral selector in the quaternary complexes. The main cleavage of $[2Phg_L+AA_{L/D}-H+Cu^{II}]^+$ involved the competitive formal Phg and $[Phg-H]^+$ losses yielding the $[Cu^{II},(Phg,AA_{L/D}-H)]^+$ and $[Cu^I,(Phg_L,AA_{L/D})]^+$ ions, respectively. In fact, the $[Phg-H]^+$ loss is not direct and requires a stepwise pathway through the prompt CO_2 release concomitant with Cu^{II} reduction. This process was therefore at the origin of the observed enantiomeric reduction. Its extent depends on the one hand, the respective chirality of the partners of the quaternary complex, and on the other hand, the basicity and acidity in the gas phase of the amino acid of unknown configuration. Concerning the thermochemical correlation, singular points appear on the one hand, for very basic and very acidic amino acids (R, K, H) and on the other hand, for D and E which are weakly basic and strongly acidic. Indeed, they exhibit a strong degradation of the enantioselective reduction.

Conclusions

The use of particular reference such as Phg allows to yield distinction of the chirality of $AA_{L,D}$ by two means either from chiral effects on the diastereomeric quaternary coppered complex dissociations described by the Cooks works or by the observed enantioselective reduction. The lowering yields of the latter for certain AAs, was explained either by the presence of zwitterion for R,H and K or by the particular reactivity of mobilizable proton on D and E in canonic form.

Novel Aspect

AA configuration by enantioselective reduction by cleavage of quaternary complexes with Phg as reference can be explained by considering AA canonic/zwitterion structures

References

1. W.A. Tao, D. Zhang, E.N. Nikolaev, R.G. Cooks. Copper(II)-assisted enantiomeric analysis of d,l-amino acids using the kinetic method: chiral recognition and quantification in the gas phase, *J. Am. Chem. Soc.* **122**, 10598–10609 (2000)
2. C. Afonso, D. Lesage, F. Fournier, V. Mancel, J.-C. Tabet. Origin of enantioselective reduction of quaternary copper d,l amino acid complexes under vibrational activation conditions. *Int. J. Mass Spectrom.* **312**, 185-194 (2012)

MOOr.08-5 **DIRECT ENANTIOMER-SELECTIVE MASS SPECTROMETRY OF MULTI-COMPONENT CHIRAL MIXTURES**

Maurice H.M. Janssen

MassSpecpecD BV, R&D, Enschede

Keywords: chiral, enantiomeric excess, direct MS, photoionization, instrumentation

Introduction

Mass spectrometry is chirally blind, it cannot directly distinguish the two enantiomers of chiral molecules. Simultaneous, enantiomer-specific identification of chiral molecules in multi-component mixtures is extremely challenging. Many established techniques for single-component analysis fail to provide selectivity in multi-component mixtures and lack sensitivity for dilute samples.

Methods

Here we show how enantiomers may be differentiated by Mass-Selected PhotoElectron Circular Dichroism (MS-PECD) using an electron–ion coincidence imaging spectrometer [1-2]. Following an ionizing circular polarized laser pulse, ions and electrons are detected in coincidence on their respective time- and position sensitive detectors. The Mass-Selected PECD reveals that the compound is chiral.

Results

As proof of concept, vapours containing ~1% of two chiral monoterpene molecules, limonene and camphor, are irradiated by a circularly polarized femtosecond laser, resulting in multiphoton near-threshold ionization with little molecular fragmentation. Large chiral asymmetries (2–4%) are observed in the mass-tagged photoelectron angular distributions. These Mass-Selected PECD asymmetries switch sign according to the handedness (R- or S-) of the enantiomer in the mixture and scale with enantiomeric excess of the chiral molecule [3-5].

Conclusions

The results demonstrate that direct MS identification of mixtures of chiral molecules is possible using the ultra-sensitive and selective technique of laser-based MS-PECD.

Novel Aspect

Direct quantitative determination of enantiomeric excess in chiral mixtures is achieved in a table-top instrument using MS-PECD. No prior chiral chromatography or enantiomer-selective clustering is needed.

References

1. C.S. Lehmann *et al.*, *J. Chem. Phys.* **139**, 234307 (2013)
2. M.H.M. Janssen and I. Powis, *Phys. Chem. Chem. Phys.* **16**, 856 (2014)
3. M.M. Rafiee Fanood *et al.*, *Nature Communications* **6**, 7511 (2015)

4. M.H.M. Janssen. *The Analytical Scientist*, May, 44 (2016)
5. M.H.M. Janssen and I. Powis, *Current Trends in Mass Spectrometry*, May, 16 (2017)

MOr.09 - Ambient MS and new ionization Methods

Chairs: Jentaie Shiea, Pietro Traldi

Monday

MOr.09-1 **Keynote: BREATH METABOLOMICS USING SECONDARY ELECTROSPRAY IONIZATION (SESI) - HIGH RESOLUTION MASS SPECTROMETRY**

Renato Zenobi

ETH Zurich, HCI E 329 - CH-8093 Zurich, Switzerland

Keywords: exhaled breath analysis • ambient ionization • secondary electrospray ionization mass spectrometry

Introduction

On-line breath analysis is a powerful approach to obtain insights into the metabolism of a person. With ambient ionization methods, this can be achieved rapidly and completely non-invasively, opening up interesting possibilities to diagnose diseases via exhaled breath, to discern diurnal changes in the metabolism, and to monitor drug use (or drug abuse) by analyzing the chemical composition of exhaled breath.

Methods

Besides SESI-MS there are other methods in use for breath analysis, notably selected ion flow tube-mass spectrometry (SIFT-MS) and proton transfer reaction-mass spectrometry (PTR-MS). These will be compared in terms of analytical performance (mass resolution, upper mass limit, detection limit, softness of ionization, etc.) and the coverage of biomarkers (polarity, volatility, molecular weight, etc.).

Results

SESI-MS affords sub-ppt limits of detection, and analysis of compounds with molecular weights up to 1000 Da. If coupled to a high-resolution instrument, SESI-MS offers excellent coverage of the breath metabolome and deep insight into the metabolic state of a person. Moreover, the high resolution allows detailed characterization of compounds that can be used as biomarkers [1]. This keynote lecture will present several examples in medical diagnosis, including the detection of novel biomarkers for respiratory diseases such as obstructive sleep apnea (OSA) [2], chronic obstructive pulmonary disease (COPD) [3], pulmonary fibrosis, and asthma. Monitoring of drug compliance and pharmacokinetics [4] via real-time SESI-MS will also be shown.

Conclusions

A variety of biomarkers that describe diseases have been detected in exhaled breath and identified via SESI-MS, often aided by chromatography, of exhaled breath condensate and comparison against commercially available or synthetic reference compounds.

Novel Aspect

Noninvasive monitoring of the detailed chemical composition of exhaled breath via SESI-HRMS.

References

1. M.T. Gaugg, et al., *Anal. Chem.* **89**, 10329-10334 (2017).

2. E.I. Schwarz et al., *Thorax* **71**, 110-117 (2016).
3. L. Bregy et al., *Clin. Mass Spectrom.* **7**, 29-35 (2018).
4. X. Li et al., *Angew. Chem. Int. Ed.* **54**, 7815-7818 (2015).

MOr.09-2 DIRECT DETECTION OF SMALL *n*-ALKANES AT SUB-PPBV LEVEL IN BREATH AIR BY PHOTOELECTRON INDUCED O₂⁺ CATION CHEMICAL IONIZATION MASS SPECTROMETRY AT kPa PRESSURE

Yan Wang - Lei Hua - Haiyang Li

Dalian Institute of Chemical Physics, Chinese Academy of Sciences, CAS Key Laboratory of Separation Sciences for Analytical Chemistry, Dalian

Keywords: small *n*-alkanes, high-pressure photoionization induced O₂⁺ cation chemical ionization, kPa pressure, direct detection, time-of-flight mass spectrometry

Introduction

Alkanes are important target analytes in the fields of petrochemical industry, environmental detection, as well as being a significant product in human metabolic processes.[1] However, online and real-time measurement of saturated hydrocarbons, especially small *n*-alkanes, using direct mass spectrometry, remains a great challenge due to their low basicity and fewer ionizable functional groups.

Methods

A high-pressure photoionization induced O₂⁺ cation chemical ionization (HPPI-OCI) source was developed and coupled with a time-of-flight mass spectrometer. An efficient online pretreatment system was developed for eliminating the interference of humidity and carbonyl compounds. High-intensity O₂⁺ reactant ion could be generated by photoelectron ionization for efficient ionization of *n*-alkanes.[2]

Results

The HPPI-OCI mass spectrum pattern of C3-C6 *n*-alkanes differed a lot at different ion source pressure. The quasi-molecular ions [M-H]⁺ were gradually dominating and the fragmentation of *n*-alkanes by O₂⁺-CI was depressed while the ion source pressure was elevated from 88 to 1080 Pa, with more than three orders of magnitude improvement in signal intensity. This was attributed to collisional deactivation effect, by which [M-H]⁺ tend to be more stable other than fragmentation, as the excess energy in the quasi-molecule ions could be more efficiently quenched at higher pressures. The achieved LODs (S/N=3) were down to 0.07-0.14 ppbv for C3-C6 *n*-alkanes in one minute. The HPPI-OCI TOFMS was successfully used in the analysis of exhaled small *n*-alkanes of healthy smokers and nonsmokers as well as the concentration variation of exhaled small *n*-alkanes after alcohol consumption.

Conclusions

This work not only provides new insights for controlling the chemical ionization process by adjusting the ion source pressures but also prompts an innovative and high-performance direct MS technique for small *n*-alkanes analysis. It is likely that the HPPI-OCI MS will be a powerful tool in various research fields such as on-site environmental analysis and high-throughput medical diagnostics.

Novel Aspect

Small *n*-alkanes can be soft ionized with O₂⁺ cation at kPa pressure and direct detected at sub-ppbv level.

References

1. Amann A., Miekisch W., Schubert J., Buszewski B., Ligor T., Jezierski T., Pleil J., Risby T., *Annual Review of Analytical Chemistry*, 7, 455-482 (2014).
2. Wang Y., Hua L., Li Q., Jiang J., Hou K., Wu C., Li H., *Analytical Chemistry*, submitted (2018).

MOr.09-3 **DIRECT COUPLING OF MAGNETIC-PARTICLE-BASED SAMPLE PREPARATION STRATEGIES TO MASS SPECTROMETRY VIA MICROFLUIDIC OPEN INTERFACE**

Marcos Tascon - Varoon Singh - Mohammad Huq - Janusz Pawliszyn

University of Waterloo, Chemistry, Waterloo

Keywords: Dispersive droplet extraction, Fast analysis, Magnetic particles, Microfluidic Open Interface, SPME-MS

Introduction

Magnetic micro- and nanoparticles are an emerging technology in the field of solid-phase extractions.[1] However, current analysis strategies are limited to off-line desorption where just a fraction of the solution is detected. Hence we developed a new design of the microfluidic open interface (MOI) which acts as an interface for the direct coupling of the magnetic-particles-based sample preparation to mass spectrometry resulting in increased sensitivity while reducing the time of analysis and sample size.[2]

Methods

Dispersive extractions of six model small molecules from droplets of PBS and urine using C₁₈ magnetic nanoparticles in droplets for 30 s were performed. On the other side, IMAC magnetic microparticles were dispersed in droplets of β -casein digests and in-droplet β -casein trypsin digestions for extraction of phosphopeptides. In both approaches, the nanoparticles were collected by applying a magnetic field followed by direct desorption via MOI and MS detection.

Results

The methodology developed for the direct coupling of magnetic particles to MS via MOI was evaluated using nano and microparticles. The desorption consisted in the generation of an isolated droplet (\approx five μ L) that wets the particles, these are held by the magnetic field applied by a 3D-printed holder. Then, the enriched droplet is rapidly drained to the ESI source where the compounds are ionized. The dispersive extraction using nanoparticles from droplets (20-40 μ L) was completed after 10 seconds. Limits of quantitation in the range 1-5 ng/mL were achieved in PBS and urine for all the compounds with an overall analysis time below 2 min. On the other hand, IMAC microparticles allowed efficient phosphopeptide extraction in less than 5 minutes from protein digests and in-droplet digest as a screening tool. Because of the high sensitivity achieved in both cases, the amount of particles used was within 30-100 μ g per sample, which is about ten times less than conventional exhaustive extraction methods.

Conclusions

A technology to efficiently couple magnetic-particle-based extractions to MS has been developed. The combination of fast desorption kinetics, low desorption volume, and no further dilutions made MOI-MS a rapid and very sensitive strategy. To proof the concept, droplet extractions using about 50 μ g of particles were tested for small and medium-size molecules. The results demonstrated good sensitivity using the minimum amount of sample and particles in a short analysis time.

Novel Aspect

A novel interface to directly and efficiently couple magnetic-particle based extractions to MS via

a microfluidic open interface was developed for fast and rapid analysis.

References

1. Vasconcelos I., Fernandes C., *Trend Anal.Chem.*, 89, 41 (2017)
2. Tascon M., Alam M.N., Gómez-Ríos G.A., Pawliszyn J., *Anal. Chem.*, 90, 2631 (2018)

MO.09-4 **DEVELOPMENT OF THE FEEDBACK CONTROL SYSTEM FOR THE ADVANCED TAPPING MODE SCANNING PROBE ELECTROSPRAY IONIZATION (T-SPESI)**

Yoichi Otsuka⁽¹⁾ - ***Bui Kamihoriuchi***⁽¹⁾ - ***Aya Takeuchi***⁽¹⁾ - ***Futoshi Iwata***⁽²⁾ - ***Takuya Matsumoto***⁽¹⁾

⁽¹⁾ *Osaka University, Department of chemistry, Toyonaka* - ⁽²⁾*Shizuoka University, Department of engineering, Hamamatsu*

Keywords: ambient sampling ionization, feedback control, scanning probe, multimodal imaging, capillary

Introduction

Ambient sampling and ionization (ASI) technique has been attracted attention due to its capability of rapid chemical analysis and/or chemical imaging. We have developed the technique named as tapping-mode scanning probe electrospray ionization (t-SPESI) [1-4], in which the unique fusion of the key techniques used in atomic force microscopy (AFM) and mass spectrometry (MS) is achieved, to realize ASI with a single capillary probe.

Methods

For the advanced SPESI, we developed the feedback control system. The optical lever method is based on the technology used in AFM. The laser light is irradiated to the capillary probe and the displacement of the reflected light is detected with the photodiode to obtain the oscillation amplitude of the probe. The position of Z-stage is controlled to maintain the oscillation amplitude during ASI with proportional-integral-derivative (PID) algorithm.

Results

The optical lever method developed in this study enables us to measure the oscillation amplitude of capillary probe. The analog to digital converter, field programmable gate array, real-time OS and PC are controlled by the homemade software coded in LabVIEW. We confirmed the oscillation amplitude is proportional to the input voltage to the piezo actuator which is attached to the capillary probe for the stable oscillation at resonant frequency. The feedback control system maintains the oscillation amplitude at the desired set point by adjusting the Z position of sample during the ASI by t-SPESI. This system can be utilized not only to stabilize the probe scanning on the sample surface with micro asperity but also to measure the surface roughness. The measurements of the mass spectra and the output signal to the z-stage is executed simultaneously. The ion images and the topographic image are reconstructed.

Conclusions

The information of the probe oscillation is utilized to improve the reproducibility of ASI with SPESI. We demonstrated the simultaneous visualization of the chemical distribution and the surface shape. This simple optical technique is expected to provide the mechanical properties of the nano-litter-volume liquid bridge containing extracted molecules from the sample.

Novel Aspect

The stability improvement of ASI and the multimodal imaging are realized by the optical lever method with t-SPESI.

References

1. Otsuka Y., Shide S., Naito J., Kyogaku M., Hashimoto H., Arakawa R., *Rapid Commun. Mass Spectrom.* **26**, 2725 (2012).
2. Otsuka Y., Naito J., Kyogaku M., Hashimoto H., Arakawa R., *Analyst* **139**, 2336 (2014).
3. Otsuka Y., Satoh S., Naito J., Kyogaku M., Hashimoto H., *J. Mass Spectrom.* **50**, 1157 (2015).
4. Kohigashi T., Otsuka Y., Shimazu R., Matsumoto T., Iwata F., Kawasaki H., Arakawa R., *Mass Spectrometry (Tokyo)* **5**, S0054 (2016).

MOOr.09-5 DIRECT MASS SPECTROMETRY ANALYSIS USING IN-CAPILLARY IONIC LIQUIDS-BASED DISPERSIVE LIQUID-LIQUID MICROEXTRACTION AND SONIC-SPRAY IONIZATION

Qiang Ma

Chinese Academy of Inspection and Quarantine, Beijing

Keywords: In-Capillary Ionic Liquids-Based Dispersive Liquid-Liquid Microextraction; Orbitrap Mass Spectrometry; Sonic-Spray Ionization; Direct Mass Spectrometry Analysis

Introduction

Public concern on the environmental and human impact of chemical approaches has encouraged analytical chemists to introduce safer and cleaner alternatives. Ionic liquids (ILs) are showing increasingly promising perspectives as green solvents in a diversity of fields. The unique physicochemical properties make ILs favorable alternative media to replace organic solvents for extraction, separation and detection techniques. In this study, a novel ILs-based dispersive liquid-liquid microextraction (DLLME) method coupled with sonic-spray ionization (SSI) has been developed for direct mass spectrometry analysis of environmental water samples. An in-situ formed dicationic IL was introduced to not only act as an effective extraction medium, but also create a positively charged complex with target analytes, allowing a significant increase in sensitivity.

Methods

A disposable borosilicate glass capillary (1.5 mm o.d., 0.86 mm i.d., 8 cm length) with a pulled tip (5 μm i.d.) was prepared using a P-1000 micropipette puller (Warner Instrument, Hamden, CT, USA). Four liquid plugs were sequentially injected into the capillary from the open end: 2 μL of aqueous $\text{Li}[\text{NTf}_2]$ solution at 30 μM , 20 μL of environmental water sample, 2 μL of aqueous $[\text{C}_4(\text{MIM})_2]\text{Br}_2$ solution at 20 μM , and 10 μL of dichloromethane. An in-capillary liquid-liquid microextraction was facilitated by the simple action of tilting the capillary up and down. The dichloromethane phase migrated to the tip end of the capillary due to higher density. The process does not need high voltage and was performed at ambient temperature.

Results

Perfluorinated compounds (PFCs), an exemplary class of chemical substances of very high concern for regulatory control due to their persistent, bioaccumulative and toxic properties, were selected for method development and validation. The environmental water plug was sandwiched between the two miscible, hydrophilic IL plugs of $\text{Li}[\text{NTf}_2]$ and $[\text{C}_4(\text{MIM})_2]\text{Br}_2$. Within the combined aqueous phase, $\text{Li}[\text{NTf}_2]$ acted as an ion-exchange reagent and underwent an in situ metathesis reaction with $[\text{C}_4(\text{MIM})_2]\text{Br}_2$ to form a new hydrophobic IL of $[\text{C}_4(\text{MIM})_2][\text{NTf}_2]_2$. The exchange process generated a cloudy, homogeneous solution with fine microdroplets of the resulting extractant phase which greatly increased the contact surface area of the IL extraction solvent resulting in enhanced extraction efficiencies for perfluorooctanoic acid (PFOA) and perfluorooctanesulfonic acid (PFOS). After hand swing, the 10 μL dichloromethane liquid slug migrated to the tip end of capillary, extracting the hydrophobic compounds. The Venturi effect occurs when the high-velocity nebulizing gas flows through the capillary tip, causing a reduction

in fluid pressure and self-pumping effect. The signal profile proved to be stable monitored by means of base peak chromatogram (BPC) and extracted ion chromatogram (EIC). The dichloromethane phase gives a clear mass spectrum with only two characteristic ions for m/z 719.1 (combination of $[C^{2+}+PFOS]^{+}$) and m/z 500.1 (combination of $[C^{2+}+Tf_2N]^{+}$). This method detect the positively charged adduct which composed of dications and analyte anions. Especially when determining the concentration of PFCs which can hardly undergo fragmentation under negative MS/MS mode. The adduct lost either a proton or a methyl imidazolium group of dications, resulting in a singly charged fragment for detection. Quantitative analysis was performed with internal standard calibration curve method, and the limits of detection (LODs) and limits of quantitation (LOQs) were 0.03 and 0.15 ng/mL for PFOS, and 0.05 and 0.21 ng/mL for PFOA.

Novel Aspect

An ILS-based dispersive liquid-liquid microextraction (DLLME) method coupled with sonic-spray ionization (SSI) has been developed for direct mass spectrometry analysis.

References

1. M. J. Trujillo-Rodriguez, V. Pino, A. M. Afonso, Trends Anal. Chem, 51, 87-106 (2013)
2. Z. Takáts, J. M. Wiseman, R. G. Cooks, Science, 306, 471-473 (2004)
3. Sebastiaan F. Teunissen, Marcos N. Eberlin, Rosana M. Alberici, Trends Anal. Chem, 90, 135-141 (2017)

MOr.10 - Polymers

Chairs: Marek M. Kowalczyk, Leopoldo Ceraulo

MOr.10-1 **Keynote: MASS SPECTROMETRY METHODS FOR THE CHARACTERIZATION OF FUNCTIONAL MATERIALS**

Chrys Wesdemiotis

Departments of Chemistry and Polymer Science, The University of Akron, Akron, OH 44325, USA

Keywords: Synthetic polymers; polymer sequence; polymer architecture; branched polymers; surface analysis

Introduction

Functional materials comprise synthetic macromolecules with well-defined structures, optimized to impart properties needed for specific biomedical, technological, or environmental applications. Multidimensional mass spectrometry (MS) methods can provide accurate and sensitive information about the molecular-level composition, structure, and architecture of such materials, a requirement for understanding and optimizing their applications.

Methods

MALDI- and ESI-MS interfaced with tandem mass spectrometry (MS²), liquid chromatography (LC) fractionation, and/or ion mobility (IM) separation were utilized to elucidate sequence and architecture. Mild thermal degradation with an atmospheric solids analysis probe (ASAP) was used to analyze crosslinked materials. Conversely, surface-layer (SL) MALDI-MS was employed to characterize and image the molecular composition and defects of solid film surfaces.

Results

MS² via collision-activated dissociation and/or electron transfer dissociation has been successfully used to determine the sequence of poly(N-isopropyl amide-methacrylic acid) copolymers, read the sequence of monodisperse copolyesters, and elucidate the architecture of

multi-arm polyethers. For more complex systems, such as branched polyglycerol conjugates, adding LC and IM separation was essential for deciphering branching architecture and degree of conjugation. Crosslinked hydrogels required mild thermal degradation via ASAP-MS; this method mainly breaks weak junctions in the polymer network leading to signature ions for both the crosslinks as well as defects such as unreacted ends. On the other hand, molecular information on the surface composition of solid films (<2 nm depth) could be gained by SL-MALDI-MS; this solvent-free technique allowed to image surface defects from foreign materials, material absence, mechanical scribing, and solvent perturbation on thin films made of poly(methyl methacrylate) and polystyrene.

Conclusions

Multidimensional mass spectrometry methods encompassing different combinations of mass analysis, tandem mass spectrometry, and orthogonal separation by polarity and/or shape enable sensitive and precise analysis of the microstructures of synthetic functional materials that are widely used in industrial, biomedical, and technological applications, but challenging to characterize by other analytical techniques due to their complexity, size, or topology.

Novel Aspect

Complete molecular-level analysis of functional materials for meaningful structure/property correlations and continuous improvement of their macroscopic properties.

MOOr.10-2 UNRAVELLING THE ARCHITECTURE OF BRANCHED POLYESTER OLIGOMERS BY COMBINED TRAPPED ION-MOBILITY SPECTROMETRY – TANDEM MASS SPECTROMETRY

Robert Voeten⁽¹⁾ - **Put van de Bram**^(2,3) - **Jan Jordens**⁽²⁾ - **Ynze Mengerink**⁽²⁾ - **Ron Peters**^(4,5) - **Rob Haselberg**⁽¹⁾ - **Peter Schoenmakers**⁽⁵⁾ - **Govert W. Somsen**⁽¹⁾

⁽¹⁾ AIMMS Division of BioAnalytical Chemistry, Department of Chemistry and Pharmaceutical Sciences, Vrije Universiteit Amsterdam, de Boelelaan 1085, 1081 HV Amsterdam, The Netherlands - ⁽²⁾ DSM Materials Science Center, Urmonderbaan 22, 6167, MD Geleen, The Netherlands - ⁽³⁾ TI-COAST, Science Park 904, 1098 XH Amsterdam, The Netherlands - ⁽⁴⁾ DSM Coating Resins, Sluisweg 12, 5145 PE Waalwijk, The Netherlands - ⁽⁵⁾ HIMS-Analytical Chemistry group, University of Amsterdam, Science Park 904, 1098 XH Amsterdam, The Netherlands

Keywords: Ion Mobility Spectrometry; Branched Polymers

Introduction

Synthetic polymers have many uses due to their large variety of potential properties. As alterations of the molecular structure, e.g. branching, result in significantly different solid-state properties, it is essential to characterize the architecture of produced polymers. We demonstrate trapped ion-mobility spectrometry with quadrupole time-of-flight mass spectrometry (TIMS-QTOFMS) for structural elucidation of branched isomeric polyester oligomers.

Methods

Polyesters of terephthalic acid (TPA) and propylene glycol (PG) with trimethylolpropane (TMP; branching agent) were synthesized. Polymer solutions (20 ppm) in THF were analysed using a Bruker timsTOF instrument allowing recording of mobilograms while continuously gathering high-resolution MS/MS data. Supercharging was realized by adding *m*-nitrobenzyl alcohol. ¹H/¹³C-NMR was used to verify polymer branching.

Results

Analysis under optimized TIMS conditions resulted in observation of multiple conformers per polymer mass. Mobility separation was achieved as branched polymers exhibit a reduced rotationally averaged 3D-projection. Tandem MS of mobility-separated isomeric polymer species

allowed their structural assignment. However, more peaks than possible structures were observed implying influence of the adduct-ion location on the respective mobilities. In verification of the latter, localization of charge by a trivalent ion, i.e. lanthanum, was achieved. For comparison, supercharging by *m*-nitrobenzyl alcohol was used to achieve multiple sodium adduct formation. Supercharging also allowed simultaneous low and high mass monitoring, showing potential towards increasing separation power. The latter provided a more accurate evaluation of the peak areas in the mobilograms, providing a semi-quantitative tool for assessment of the branching content of the polymer samples.

Conclusions

The complex and diverse properties of industrially relevant polyesters provided an ideal sample to show the added value of TIMS for polymer characterization. Structural elucidation of the multitude of isomers is realized by combining mobility traces and MS/MS data. Supercharging allowed assessing the branching content semi-quantitatively, and enabled mass observation of three times the number average molecular mass.

Novel Aspect

Use of TIMS on supercharged polyesters allows assessment of the branching content, structure and adduct location of polyester oligomers.

Acknowledgement

This research was carried out within the HOSAna project which is funded by the Netherlands Organization for Scientific Research (NWO) in the framework of the Programmatic Technology Area PTA-COAST4 and is supported by Bruker, DSM and Postnova Analytics (Project number 053.21.117).

MOr.10-3 *MS/MS READING OF POLYURETHANE BARCODES USED AS TAGGANTS IN POLYMERIC MATERIALS*

Salomé Poyer⁽¹⁾ - *Supansa Chimjarn*⁽¹⁾ - *Denise Karamessini*⁽²⁾ - *Jean-François Lutz*⁽²⁾ - *Laurence Charles*⁽¹⁾

⁽¹⁾ Institut de Chimie Radicalaire, Aix-Marseille Université, Marseille - ⁽²⁾ Institut Charles Sadron, Université de Strasbourg, Strasbourg

Keywords: Sequence-defined polymers, tandem mass spectrometry, molecular barcodes

Introduction

Advanced polymerization techniques allow production of synthetic chains with a perfect control of the molecular structure and of the comonomeric sequence [1]. One key application of such sequence-defined polymers is information storage: a message can be chemically encoded in the chain by using a set of two comonomers and can be read using tandem mass spectrometry (MS/MS)[2]. These polymers can be used as molecular barcodes in anti-counterfeiting applications [3].

Methods

Polyurethane (PU) oligomers were synthesized as previously described [2]. Reading of different sequence-defined polyurethanes was evaluated in electrospray tandem mass spectrometry from deprotonated ions using a Q-TOF instrument. Extraction procedures of PUs were optimized as a function of the nature of the membrane.

Results

Sequence-defined polyurethanes (PU) are composed of carbamate-based comonomers in which the coded moiety is a small alkyl segment that contains (1-bit) or not (0-bit) methyl substituents. These are extremely easy to read in the negative mode due to i) single cleavage of the carbamate bond in each monomer and ii) unique location of the negative charge at the α chain-end. As a result, the message coded in these PU chains can easily be read by MS/MS as the distance between two fragment ions corresponds to the mass of one or the other comonomer. Due to such a high MS/MS readability, PUs were tested as molecular barcodes for the labeling of different membranes made of PS, PVC, PET and PVA. Optimization of the PU extraction step for each polymeric substrate permitted to reach detection limits in the ppm range for MS/MS sequencing. Moreover, the tagging homogeneity of membranes was evaluated from ESI-MS experiments.

Conclusions

Sequence-defined PUs were used as molecular taggants for the labeling of various polymeric materials. Thanks to the renowned high specificity and sensitivity of MS/MS, taggants were recovered and identified even at very low concentration levels. Selectivity of taggant extraction from the membranes prevented any matrix effects during the electrospray ionization step, allowing direct sample introduction and hence fast reading of the molecular code.

Novel Aspect

PU barcodes were successfully recovered and identified by MS/MS from polymeric materials in which they were used as anti-counterfeiting labels

References

1. Lutz J.-F., Lehn J.-M., Meijer E. W., Matyjaszewski K., *Nat. Rev. Mater.* 1, 16024 (2016).
2. Amalian J.A., Poyer S., Petit B.E., Telitel S., et al., *Int. J. Mass Spectrom.* 421, 271-278 (2017).
3. Karamessini D., Poyer S., Charles L., Lutz J.-F., *Macromol. Rapid Commun.* 38, 1700426 (2017).

MOr.10-4 ELECTROSPRAY IONIZATION MASS SPECTROMETRY OF NONPOLAR POLYISOBUTYLENES

Sándor Kéki⁽¹⁾ - **Lajos Nagy**⁽¹⁾ - **Tibor Nagy**⁽¹⁾ - **György Deák**⁽¹⁾ - **Ákos Kuki**⁽¹⁾ - **Mihály Purgel**⁽²⁾ - **Mijid Narmandakh**⁽³⁾ - **Béla Iván**⁽⁴⁾ - **Miklós Zsuga**⁽¹⁾

⁽¹⁾ University of Debrecen, Department of Applied Chemistry, Debrecen – ⁽²⁾ University of Debrecen, Department of Applied Chemistry, Debrecen – ⁽³⁾ University of Debrecen, Department of Physical Chemistry, Debrecen – ⁽⁴⁾ Research Center for Natural Sciences, Hungarian Academy of Sciences, Polymer Chemistry Research Group, Institute of Materials and Environmental Chemistry, Budapest – ⁽⁵⁾ Research Center for Natural Sciences, Hungarian Academy of Sciences, Research Group, Institute of Materials and Environmental Chemistry, Budapest

Keywords: ESI-MS, polyisobutylenes, adduct ions, broadband MS/MS, survival yield

Introduction

Polyisobutylenes (PIBs) is an important class of non-polar hydrocarbon polymers produced by cationic polymerization of isobutylene [1]. Owing to their outstanding chemical resistance and biocompatibility, PIBs can be applied in many fields spanning from coatings to various biomaterials [1]. In this lecture, we present a method based on the adduct ion formation of PIBs with anions [2], called anion-attachment method [3] for the characterization of low molecular weight PIBs using electrospray ionization mass spectrometry (ESI-MS).

Methods

Polyisobutylenes with different molecular weights and end-groups were synthesized by living cationic polymerization. Measurements were performed with a MicroTOF-Q type Qq-TOF MS

instrument (Bruker Daltonics, Bremen, Germany). The critical (activation) energies for the dissociation of $[\text{PIB}+\text{anion}]^-$ into PIB + anion were determined by means of broadband MS/MS measurements using the survival yield (SY) method and supported by high-level quantum chemical calculations (DFT).

Results

The $[\text{PIB}+\text{anion}]^-$ adducts were generated using electrospray ionization (ESI) by adding the corresponding salts to the PIB solution. It was found that independently of the end-groups, presence or absence of aromatic initiator moiety, the PIB derivatives investigated in this study were capable of forming adduct ions with various anions such as NO_3^- and Cl^- ions, thus allowing the characterization of these compounds using the negative ion mode of ESI-MS. The dissociation of the $[\text{PIB}+\text{anion}]^-$ adduct was shown to yield neutral PIB and charged anion and that the collision energy necessary to obtain 50 % dissociation of the precursor ion (CE_{50}) changed linearly with the number of isobutylene units. Furthermore, it was found that the binding enthalpies of the halide ions to PIB, determined by using broadband MS/MS and the SY methods, increased in the order of $\text{I}^- < \text{Br}^- < \text{Cl}^- < \text{F}^-$. This order is in line with that obtained by DFT calculations and correlate with the size of the halide ion.

Conclusions

The anion-attachment under ESI conditions proved to be an appropriate method for ionizing the non-polar polyisobutylene derivatives. Furthermore, applying broadband MS/MS and the SY method information on the strength of the interaction between the respective anion and the non-polar polymer chain can be gained.

Novel Aspect

Our work has revealed that polymers even with high hydrophobicity can be effectively characterized by ESI-MS which is known to be mainly capable of investigations of polar compounds.

Acknowledgement

The work was supported by the GINOP-2.3.2-15-2016-00041 project and by the grant K-116465.

References

1. Varion J. P., Spassky N.: *Industrial Cationic Polymerizations* Marcel Dekker, New York, Basel (1996)
2. Nagy L., Nagy T., Deák Gy., Kuki Á., Purgel M., Narmandakh M., Iván B., Zsuga M. Kéki S. J. *Am. Soc. Mass Spectrom.* 27(3), 432-442 (2016)
3. Nagy L., Kuki Á., Deák Gy., Purgel M., Vékony Á., Zsuga M., Kéki S. J. *Phys. Chem. B.*, 120, 9195-9203 (2016)

MOr.10-5 ANALYSIS OF POLYSTYRENE BASED MICROPLASTICS IN THE ENVIRONMENT

Gabriella Francesca Schirinzi - Marta Llorca - Marinella Farré - Damià Barceló

CSIC, IDAEA, Barcelona

Keywords: polystyrene, microplastics, seawater, HPLC-APPI/HRMS

Introduction

Marine anthropogenic litter is a severe environmental problem[1]. One of the main issues is the extreme stability of plastic wastes[2], which can generate nanoplastics (NPLs) and microplastics (MPLs)[3]. The analysis of these plastic particles should consider a combination of different

analytical tools for the quantitative/qualitative analysis of MPLs/NPLs. These studies have been carried out using polystyrene (PS) as a representative polymer.

Methods

A standard of PS (MW ~1.2kDa) was solubilized in toluene at room temperature. For the mass spectrometric analysis, the performance of different ionisation sources (ESI, APCI, APPI, DART, MALDI) were optimized and assessed under negative and positive conditions. The chromatographic separation was achieved by size-exclusion with a GPC-column and using toluene as mobile phase. Mass spectrometric analyses were carried out using a QExactive mass spectrometer.

Results

The optimal results were obtained by HPLC-HRMS equipped with atmospheric pressure photoionization (APPI) source operating in negative mode. The ionisation of PS resulted in a characteristic profile of signals, each of them separated by 104 m/z. The addition of a chromatographic separation with GPC column showed a good performance for qualitative and semi quantitative measurements in complex extracts. The APPI source was by far the most efficient ionisation source for this polymer in comparison with ESI and APCI.

Additionally, the potential of ambient air mass spectrometry techniques, such as DART, for the detection and identification of PS was assessed. The results suggested that MALDI and ambient air ionization techniques result in low ionisation efficiencies in comparison with APPI.

These results were combine with others different techniques such as TGA and FT-IR which are useful tools for the characterisation and identification of polymers but not to quantify them.

Conclusions

GPC-LC coupled to APPI-QExactive is a suitable approach for the analysis of low molecular weight PS (< 4000Da) and applies to real seawater samples. However, the use of LC-APPI-HRMS should be complemented by other techniques such as TGA, DSC and FT-IR, which allow obtaining qualitative and semi-quantitative information about high-molecular weight polymers, which may also be present in the environment.

Novel Aspect

The environmental analysis of polymers should rely on a protocol, which combines different types of techniques for the complete characterization of polymers.

References

1. Bergmann, M., L. Gutow, and M. Klages, *Marine Anthropogenic Litter*. 2015: Springer.
2. Li, W., H. Tse, and L. Fok, *Plastic waste in the marine environment: A review of sources, occurrence and effects*. *Science of The Total Environment*, 2016. **566**: p. 333-349.
3. Carr, S.A., J. Liu, and A.G. Tesoro, *Transport and fate of microplastic particles in wastewater treatment plants*. *Water research*, 2016. **91**: p. 174-182.

TUESDAY, August 28th, 2018

PL03: Plenary Lecture: Thomson Medalist

Chair: Catherine E. Costello

Albert J. R. Heck

Utrecht University, The Netherlands

Keywords: native mass spectrometry, hybrid fragmentation techniques, mass spectrometry based structural biology, circadian clock, ribosomes

Introduction

Around for more than a century mass spectrometry (MS) is blooming more than ever and applied in all aspects of the natural and life sciences. In the last two decades MS has become routine for the high-throughput analysis of proteins and their interactions and post-translational modifications. A next phase will be to gather in-depth structural and functional information on all these proteins and their cellular assemblies, requiring novel MS-based approaches.

Methods

To appreciate protein function requires structural models, and detailed information on how these evolve during cellular perturbations. MS can also contribute to this next level of protein analysis, giving rise to an emerging new area in structural proteomics. Our group is constantly developing methods using native MS, HD-exchange MS and cross-linking MS with the aim to bridge the gap between interaction proteomics and structural biology.

Results

In the first part of the talk novel developments in MS instrumentation for native MS will be described, especially new modifications to an Orbitrap based instrument that offers high-sensitivity and mass resolution, allowing an in-depth detailed analysis of highly glycosylated glycoproteins, protein assemblies up to even whole intact ribosomes and viruses [1]. We present high-resolution native MS data on ribosomal particles from *E. Coli*, spinach and human cells, and by combining this data with top-down and bottom-up proteomics analysis come to an unprecedented detailed structural insight into the micro-heterogeneity observed in these particles [2].

The second part of the talk will focus on an intriguing protein assembly involved in the circadian timing in cyanobacteria. Herein, by using a combination of native, HD exchange MS, XL-MS and cryoEM, we could define a novel structural model improving our understanding of the evolutionary oldest system controlling circadian rhythms [3].

Conclusions

Through advances in technologies the role of MS in structural and molecular biology is rapidly expanding [4]. Especially hybrid MS approaches combining top-down and bottom-up proteomics analysis with native MS, HD exchange MS or XL-MS will assist us to provide unprecedented detail on the structural and function heterogeneity present in dynamically evolving protein assemblies such as intact ribosomes and viruses.

Novel Aspect

The here described hybrid MS approaches will become an integral part of structural biology, complementing cryo-EM and molecular modelling.

References

1. van de Waterbeemd M., Fort K.L., Boll D., Reinhardt-Szyba M., Routh A., Makarov A., Heck A.J.R., *Nat Methods*. 14, 283-286 (2017)
2. van de Waterbeemd M., Tamara S., Fort K.L., Damoc E., Franc V., Bieri P., Itten M., Makarov A., Ban N., Heck A.J.R., *Nat Commun*. 9, 2493 (2018)

TOr.01 - Therapeutic drug monitoring and drug discovery

Chairs: Gérard Hopfgartner, Fabio Garofolo

TOr.01-1 **Keynote: THERAPEUTIC DRUG MONITORING AND DRUG DISCOVERY: IS IT ALL KNOWN OR ARE WE INNOVATIVE AND GAME CHANGING?**

Ron M. A. Heeren⁽¹⁾ - **Filip Cuyckens**⁽¹⁾ - **Rob J. Vreeken**^(1,2)

⁽¹⁾ Maastricht MultiModal Molecular Imaging (M4I) institute, Maastricht, The Netherlands - ⁽²⁾ Janssen R&D, Beerse, Belgium

Keywords: Quantitation, omics, instrumentation, drug discovery

Introduction

Nowadays, mass spectrometry, mostly hyphenated with either GC or LC, is THE prime technology applied in therapeutic drug monitoring (TDM) and drug discovery (DD). The accurate quantitation of drugs and related metabolites are essential for personalized treatment scenario's in the clinic, while both quantitative and qualitative data of these drugs and associated metabolites or PD effects are required in drug discovery.

Methods

In order to address all questions in both TDM and DD a wide range of instruments, spanning from single quads till high resolution MS-systems, like e.g. the Orbitrap and even FTICR systems, are applied in conjunction with state-of-the-art separation sciences to provide reliable data supporting decision making in either continuation or stopping further development or treatment.

Results

Here we address several of the challenges we are facing in the two application areas and present some of the solutions. Next to small molecules, more and more larger molecules, both as proteins but also Post-Translational Modifications (PTM) need to be monitored, quantified and related to either beneficial or adverse effects from initial PK until lead optimization and further development. Also in TDM, next to the traditional small molecules, proteins, antibodies, drug-antibody complexes etc. are getting more attention in the clinic. All these new area's are relying on sensitivity, specificity and both spectral and chromatographic resolution. Hence, a diversity of applications will be addressed. Besides, the growing interest in applying not only targeted but also omics-like strategies to profile off-target effects, pharmacodynamics and drug related compounds is gaining popularity in early screening already. In TDM on the other hand, these technologies are far from routine at the moment.

Conclusions

Although, new technology, either being hardware, software is developed, there still is a need for higher sensitivity, higher acquisition speed, higher dynamic range and higher resolution in both TDM and DD. Here some examples are shown addressing some of the challenges.

Next, we show some of the novel applications rising at the horizon to further improve reliable and quantitative as well as qualitative information, required to make better and faster decisions on medication or progression of novel drugs in further development.

TOr.01-2 A HIGHLY SELECTIVE AND SENSITIVE LC-MS/HRMS ASSAY FOR COPROPORPHYRIN-I AND -III, EMERGING ENDOGENOUS BIOMARKERS OF OATP, IN FIRST-IN-HUMAN CLINICAL TRIALS

Ragu Ramanathan - Amanda King-Ahmad - Christopher L. Holliman, A. David Rodrigues

Pfizer Inc., Medicine Design, Groton, CT 06340, USA

Keywords: High Resolution Mass Spectrometry (HRMS), Transporter Biomarkers, Organic Anion Transporting Polypeptides (OATP), Coproporphyrins, Drug-Drug-Interactions (DDI)

Introduction

The evaluation of transporter biomarkers in early clinical trials could help identify NCE candidates with true drug-drug interaction (DDI) liabilities and thus reduce the number of costly clinical DDI studies being conducted. Coproporphyrin-I and Coproporphyrin-III are two such endogenous isomeric compounds which have emerged as potentially viable biomarkers of organic anion transporting polypeptide (OATP) mediated DDIs.

Methods

A UHPLC-MS/HRMS assay was developed using, previously unreported, CP-I/CP-III doubly charged ions as precursor ions to improve the assay sensitivity over the singly charged precursor ions. The TOF-MS/MS selectivity of the assay was improved by monitoring and post-acquisition summing the responses from the major doubly charged product ions. The LC-MS/HRMS assay, developed and qualified, supported the quantitative determination of CP-I/CP-III in human plasma.

Results

In this study, we used both singly- ($[M+H]^+$) and doubly-charged ($[M+2H]^{2+}$) precursor ions of CP-I/CP-III respectively at m/z 655.3 and 328.1 to evaluate high-resolution MS (TOF-MS) for untargeted and MS/HRMS for targeted quantification data on a high-speed Triple TOF 6600 LC-MS/MS system. For the first time, we were able to demonstrate that the doubly-charged ($[M+2H]^{2+}$) precursor ions of CP-I/CP-III are present and detectable with 5-10 times more intensity.

Calibration and quality control samples were prepared to cover a range of 50-1000 pg/mL in 3-times charcoal stripped plasma, and parallelism to native plasma established. The LC-MS/HRMS assay, developed and qualified, supported the quantitative determination of CP-I/CP-III in human plasma with an assay range of 0.05-10 ng/mL. In control pooled human plasma, levels of CP-I and CP-III ranged from 0.45 and 1.1 ng/mL and from 50 and 500 pg/mL, respectively. In control pooled human urine, levels of CPs were at least 10-fold higher than those detected in plasma.

Conclusions

The NCE in this study had been shown to inhibit OATP *in vitro* and the probe drug DDI study was incorporated into the clinical study. At the highest dose of the NCE, CP-I showed a 1.37-fold increase in plasma AUC_{0-48hr} versus placebo, which was comparable to the plasma AUC_{0-48hr} ratio change observed with atorvastatin (1.62). These results support that UHPLC-MS/HRMS assay is sensitive and selective to measure endogenous biomarkers of DDIs.

Novel Aspect

Clinical application of UHPLC-MS/HRMS assay for endogenous OATP biomarkers, CP-I and CP-III, and correlation of AUC_{0-48hr} ratio changes with those from atorvastatin.

TOr.01-3 DEVELOPMENT OF AN EX VIVO SOLID-PHASE MICROEXTRACTION CALIBRATION MODEL TOWARDS IN VIVO QUANTITATION OF DOXORUBICIN IN LUNGS DURING A NOVEL CHEMOTHERAPEUTIC PROCEDURE.

Anna Roszkowska⁽¹⁾ - **Marcos Tascon**⁽¹⁾ - **Barbara Bojko**⁽²⁾ - **Krzysztof Goryński**⁽²⁾ - **Germán Augusto Gómez-Ríos**⁽¹⁾ - **Nikita Looby**⁽¹⁾ - **Pedro Reck dos Santos**⁽³⁾ - **Marcelo Cypel**⁽³⁾ - **Janusz Pawliszyn**⁽¹⁾

⁽¹⁾ University of Waterloo, Department of Chemistry, Waterloo - ⁽²⁾ Collegium Medicum in Bydgoszcz, Department of Pharmacodynamics and Molecular Pharmacology, Torun - ⁽³⁾ University of Toronto, Latner Thoracic Surgery Research Laboratories, Toronto General Research Institute, University Health Network and Department of Surgery, Toronto

Keywords: *in vivo* solid-phase microextraction, doxorubicin, therapeutic drug monitoring, LC-MS/MS

Introduction

The real-time monitoring of the concentration of doxorubicin (DOX) in lung tissue aims to enhance the effectiveness of applied chemotherapy while reducing severe side effects [1]. To precisely assess the level and biodistribution of DOX, *in vivo* solid-phase microextraction (SPME) has been introduced as minimally invasive technology that integrates sampling, extraction and sample clean-up, and allows for repeated sampling of different tissue compartments [2-4].

Methods

In optimized SPME conditions, sterile 15 mm (45 µm thickness) mixed-mode fibers underwent preconditioning in acetonitrile/water mixture (80:20, v/v) for 60 min. Next, the fibers were directly introduced into lung tissue for 20 min extraction of DOX under static conditions, followed by a desorption step in acetonitrile/water (80:20, v/v) with 0.1% formic acid for 60 min and finally the obtained extract was subjected to LC-MS/MS analysis.

Results

In SPME method development, several parameters that could influence the extraction efficiency of DOX (e.g., selection of optimal preconditioning and desorption solvents, fiber sterilization) were evaluated. SPME extraction in lung tissue reached equilibrium in 20 min. The entire procedure was performed without internal standard correction, enabling the methodology to be implemented in *in vivo* SPME sampling. Under optimized conditions, SPME exhibited excellent linearity ($R^2 \geq 0.99$) within the range of 2.5-50 µg/g, a LOQ of 2.5 µg/g, good precision (RSD<14.7%) and accuracy (≈ 103.2). Relatively short sampling time ensures that alterations in the level of DOX can be monitored during chemotherapy, and that a near-real-time profile of drug biodistribution in the lung tissue can be provided. Due to the minimal invasiveness of SPME probes to the living system, sampling can be performed several times in different tissue compartments during the treatment of cancer metastases in lungs with the use of *in vivo* lung perfusion (IVLP) [5].

Conclusions

An analytical system based on SPME sampling coupled to LC-MS/MS method was developed to accurately and precisely quantitate the level of the anticancer drug, DOX, by direct sampling of lungs. The optimized SPME method facilitates near to real-time monitoring of the concentration and biodistribution of DOX in the lung tissue of a living organism (*in vivo* SPME) during IVLP in the operating room.

Novel Aspect

In vivo SPME method for repeated sampling of DOX under equilibrium conditions without internal standard correction was developed to monitor drug level in different lung areas.

References

1. Tacar O., Sriamornsak P., Dass C.R., *J. Pharm. Pharmacol.* 65, 157–170 (2013).
2. Ouyang G., Vuckovic D., Pawliszyn J., *Chem. Rev.* 111, 2784–2814 (2011).
3. Bojko B., Gorynski K., Gomez-Rios G., Knaak J.M., Machuca T., Cudjoe E., Spetzler V.N., Hsin M., Cypel M., Selzner M., Liu M., Keshjavee S., Pawliszyn J., *Lab. Invest.* 94, 586–94 (2014).
4. Roszkowska A., Tascon M., Bojko B., Goryński K., Dos Santos P.R., Cypel M., Pawliszyn J., *Talanta* 183, 304–310 (2018).
5. Dos Santos P.R., Sakamoto J., Chen M., Linacre V., Arce C., Liu M., Waddell T.K., Keshavjee S., Cypel M., *Ann. Thorac. Surg.* 101, 2132–2140 (2016).

TO.r.01-4 CHARACTERIZATION OF PACLITAXEL METABOLITES BY UHPLC-MS/MS: BIOMEDICAL APPLICATION OF THE STUDY IN TARGETED DRUG DELIVERY SYSTEMS

Cristina Sottani - Elena Grignani - Emanuele Studioso - Danilo Cottica

Istituti Clinici Scientifici Maugeri, Environmental Research Center, Pavia ⁽¹⁾

Keywords: paclitaxel, metabolite characterization, UHPLC-MS/MS, nanoparticles, new drug delivery systems

Introduction

Recently has been proved that mesenchymal stem cells (MSCs) are able to package and deliver paclitaxel through their extracellular vesicles, suggesting the possibility of using MSCs to develop drugs with a higher cell-target specificity. The main objective of this study was to exploit the ability of mass spectrometry in order to perform structural studies on paclitaxel molecule because this drug follows a specific fragmentation pathway, which cleaves the molecules in two major fragment ions. These fragments are considered markers in biological fluids.

Methods

The bile was collected every hour from ½ h to 3 h. An aliquot of 300 mL was picked up from every samples, diluted with acetonitrile 1:1 (v:v) and then centrifuged. The supernatant was evaporated to dryness. The residue, dissolved in 100 mL with 0.02 M ammonium acetate buffer and acetonitrile 65:35 was directly injected into LC-MS/MS system. The entire effluent from the column of the control and pooled bile samples were analyzed in mass spectrometry. Product ion LC-MS/MS scans were obtained.

Results

It is possible to take advantage of paclitaxel mass spectrometric behavior to establish the main features of drug metabolism, given that the hydroxylation process takes place on the two major fragment ions. Paclitaxel has extensively been studied both in the laboratory and in the clinic. However, its complete metabolic pathway in rats and humans have not yet been fully understood.

We concluded that the simply observation of intense molecular ions, even if well separated from the chemical noise, is not enough to determine the presence of paclitaxel analogs as well as their chemical modifications. Tandem mass spectrometry allowed to detect the presence of 9 metabolites, distinguishing them from the other endogenous contaminants.

These metabolites were recognized as three di-hydroxytaxol, four mono-hydroxytaxol, one deacetyltaxol and one containing the taxane ring.

With this methodology, we were able to identify four new metabolites of paclitaxel belonging to the di-hydroxy and mono-hydroxy series.

Conclusions

The assessment of the most important paclitaxel metabolites represents the first structured approach to pharmacokinetic studies (PKs) of this drug used in modern cancer therapy that is based on the development of new formulations. Pharmaceutical formulations such as liposome, micelle, emulsion and nanoparticles are potentially useful for drug delivery systems. The pivotal characteristics of new drug delivery systems, such as engineered nano-materials, have to include biodegradability, biocompatibility and non-toxicity for human cells [1-2].

Novel Aspect

Paclitaxel has been recently loaded in different types of engineered nano-materials because they are demonstrated to penetrate specific tissues and reach directly the tumor site.

References

1. Danhier F., et al., PLGA-based nanoparticles: An overview of biomedical applications. *Journal of Controlled Release*, 161, 505 (2012).
2. Masood F., Polymeric nanoparticles for targeted drug delivery system for cancer therapy. *Materials Science and Engineering C*, 60, 569 (2016).

TOr.01-5 **THE EFFECTS OF LOW TEMPERATURE STORAGE ON THE STRUCTURE OF HERCEPTIN®: TO FREEZE, OR NOT TO FREEZE: THAT IS THE QUESTION**

Rosie Upton⁽¹⁾ - **Kamila Pacholarz**⁽¹⁾ - **David Firth**⁽²⁾ - **Sian Estdale**⁽²⁾ - **Perdita Barran**⁽¹⁾

⁽¹⁾ Manchester Institute of Biotechnology, University of Manchester, Manchester - ⁽²⁾ Covance Ltd., Biopharma CMC, Harrogate

Keywords: IgG1, Herceptin®, cryo-stability, IM-MS, intrinsic-fold

Introduction

For research and drug development mAbs are usually examined as fresh solutions stored between 2 – 4°C. By contrast the manufactured product is freeze-dried and stored at much lower temperatures. Here we use IM-MS and HDX-MS to probe structural differences between lots of the IgG1 mAb, Herceptin® as well as the effect of cold storage on the intrinsic conformation of this protein and its glycoforms.

Methods

DT-IM-MS measurements were used to compare refolded mAbs against the freeze-dried API. Collision induced unfolding (CIU) measurements were performed using an RF-confining linear drift field Synapt G2 along with circular dichroism (CD) to provide stability profiles. HDX-MS analyses enabled the dynamic properties to be compared. We compare data obtained from fresh, freeze-dried and reconstituted solutions.

Results

The fresh Herceptin® sample has a significantly broader collision cross section distribution (CCSD 45 – 100 nm²) compared with the frozen originator lots (60 – 80 nm²) despite being identical in terms of protein sequence. This is the case across 3 lots of Herceptin® as well as a lower purity ('research grade' trastuzumab) sample, which are extensively characterised with IM-MS as well as HDX-MS. Data from the renatured protein is compared with these two extreme cases in order to determine if freeze-storage fundamentally alters the protein fold. In each case the disulphide bridges are mapped. Global conformational stability is further probed using collisional activation prior to ion mobility analysis, and the profiles compared for lot-to-lot variability and storage effects. CD is used to examine the differences in secondary structure. As controls we have

analysed the NIST mAb standard[1] and an intact mass check standard from Waters, using the same combination of mass spectrometry and biophysical techniques.

Conclusions

We see significant differences in the conformational profile of Herceptin® depending on its treatment and between lots. This indicates that harsh storage conditions such as deep-freezing, do indeed restrict flexibility in the global conformation of Herceptin®. We present data that shows if complete denaturation of both the fresh and frozen mAb can be used to restore flexibility and consider the implications of these findings for drug processing following discovery.

Novel Aspect

Insights into the effects of cold storage upon the stability profiles, intrinsic-fold and global conformations of multiple Herceptin® lots.

References

1. "NIST Monoclonal Antibody Reference Material 8671," [Online]. Available: <https://www.nist.gov/programs-projects/nist-mono-clonal-antibody-reference-material-8671> (2018)

TOr.02 - Single cell

Chairs: Peter Nemes, Giuliana Bianco

TOr.02-1 **Keynote: SINGLE CELL MASS SPECTROMETRY: FROM QUALITATIVE TO QUANTITATIVE STUDIES**

Zhibo Yang

University of Oklahoma

Keywords: Live single cell, ambient mass spectrometry, metabolomics, microscale bioanalysis, Single-probe

Introduction

Cell is the basic unit of life. Cells are traditionally studied through population analysis, and averaged results are obtained. Understanding individual cells is critical for advances in many areas ranging from cell biology to clinical treatment. Due to its high sensitivity and broad detection range, MS is promising technique for single cell analysis.(1-4)We have developed novel techniques using microscale probe to study live single cells under ambient conditions.

Methods

We have developed a microscale sampling and ionization device (the Single-probe) that can be used for SCMS studies.(5-7)The Single-probe is fabricated by combining a laser-pulled dual-bore quartz needle (tip size smaller than a cell) with a silica capillary and a nano-ESI emitter. In SCMS experiments, cells are attached to an XYZ-translational stage for sample movement. Cell insertion is monitored using digital microscopes as the Z-stage is lifted precisely.

Results

(1) Qualitative studies. Cancer cell lines were used as model systems, and they were cultured and treated with anticancer drugs. We have successfully detected anticancer compounds, drug metabolites, and cell metabolites inside single cancer cells. (5-7)

(2) Semi-quantitative studies. Cancer stem cells (CSCs) and non-stem cancer cells (NSCCs) were used for comparison studies. Metabolomic analyses indicated that CSCs possess higher abundances of unsaturated lipids, and the level of desaturation can be altered using inhibitors.

(3) Quantitative studies. To quantify the amount of anticancer drugs in cells, sampling solvent containing the deuterated drug was delivered into a target cell to extract intracellular compounds. Both target and deuterated molecules were simultaneously detected, and abundances of target molecule were estimated from their relative ion intensities, the concentration of the internal standard, solvent flowrate, and data acquisition time. Standard HPLC/MS analysis of cell lysates were done for comparison.

Conclusions

Using the Single-probe SCMS techniques, we conducted live single cell studies including qualitative, semi-quantitative, and quantitative analyses. Broad ranges of molecules were detected. Experimental data were comprehensively analyzed to extract metabolomic information and reveal mechanisms of anticancer drugs and cell metabolism. Drug uptake amounts by single cancer cells were estimated. Our techniques can be potentially applied in clinical studies.

Novel Aspect

Multi-functional microscale device was coupled to MS for live single cell studies. Experiments and advanced data analysis were combined to obtain metabolomic information.

References

1. Miura D., Fujimura Y., Wariishi H. *Journal of proteomics*. 75, 5052-5060 (2012). Epub 2012/03/01.
2. Passarelli M. K., Ewing A. G. *Current opinion in chemical biology*. 2013. Epub 2013/08/21.
3. Lanni E. J., Rubakhin S. S., Sweedler J. V. *Journal of proteomics*. 75, 5036-5051 (2012). Epub 2012/04/14.
4. Boggio K. J., Obasuyi E., Sugino K., Nelson S. B., Agar N. Y. R., Agar J. N. *Expert Rev Proteomic*. 8, 591-604 (2011).
5. Rao W., Pan N., Yang Z. *Journal of visualized experiments, JoVE*. (112) e53911 (2016).
6. Pan N., Rao W., Standke S. J., Yang Z. *Anal Chem*. 88, 6812-6819 (2016).
7. Pan N., Rao W., Kothapalli N. R., Liu R. M., Burgett A. W. G., Yang Z. *Anal Chem*. 86, 9376-9380 (2014).

TO.r.02-2 ANALYSIS OF SINGLE MAMMALIAN CELLS WITH SCOPE-MS WITH MPOP

Harrison Specht⁽¹⁾ - **Guillaume Harmange**⁽¹⁾ - **David Perlman**⁽¹⁾ - **John Neveu**⁽²⁾ - **Christopher Adams**⁽³⁾ - **Nikolai Slavov**⁽¹⁾ - **Bogdan Budnik**⁽²⁾

⁽¹⁾ Northeastern University, Bioengineering Department, Boston - ⁽²⁾ Harvard University, FAS, Division of Science, Cambridge - ⁽³⁾ Bruker, Proteomics, San Jose

Keywords: single cell proteomics

Introduction

Quantifying proteins in single cells has a long history. For decades, scientists and physicians have used antibodies and fluorescent proteins to identify or quantify a few different proteins per cell. These methods have enabled new discoveries, clinical applications, and even spawned new fields based on measuring just a few different proteins per cell. We introduced a new technique SCOPE-MS¹ that allows measurements of hundreds of proteins in single mammalian cells. The biggest limitation of that technique was throughput of sample preparation procedure and losses

during clean up steps. To improve throughput and minimize losses we introduce a new sample preparation protocol for single cell proteomics.

Methods

A major limitation to applying quantitative mass-spectrometry to a small amount of sample, such as single mammalian cells, are the losses during sample cleanup. To relieve this limitation, we developed a Minimal Proteomic sample Preparation (mPOP)² method for culture-grown mammalian cells. mPOP obviates cleanup, and thus eliminates cleanup related losses while simplifying and expediting sample preparation for analysis of single cells.

Results

SCOPE MS technique with mPOP preparation procedures were applied for several types of single cells—analysis. As previously, new preparation procedures were tested for differentiation of Jurkat vs U cells to prove that new preparation protocols did not aggravate the ability to accurately quantify proteomes of the single cell. This new technique was applied to analysis of human activated Jurkat cells, mouse brain neurons and astrocytes. Obtained data shows that the new method was able to do quantitative analysis of hundreds of mammalian cells and show clustering of different types of cells according to their unique proteome properties.

Conclusions

Using mPOP with cell-sorting and liquid handling, we can lyse thousands of single cells in minutes, and quantify more than 700 proteins across 100 single cells in less than 12 hours. Likewise, mPOP enables protein measurements in 10, 100, and 1000 cell samples with unprecedented depth and throughput.

Novel Aspect

SCOPE MS with mPOP preparation procedure unlocks our ability to quantify hundreds of cells per day, with depth up to a thousand proteins. This type of analysis will allow researchers to create better models of protein-protein interaction and build new models of biological systems as well as new biomarkers of particular cell types.

References

1. Budnik, B., Levy, E., Harmange, G. & Slavov, N. Mass-spectrometry of single mammalian cells quantifies proteome heterogeneity during cell differentiation. *bioRxiv* 1, DOI: 10.1101/102681 (2017).
2. Specht H, Harmange, G, Budnik B, Slavov N, A minimal cell lysis and bottom-up proteomics sample preparation (mPOP) technique for mammalian cells in suspension. Manuscript in preparation 2018.

TOr.02-3 **LIVECELL MS METHODOLOGY IN THE CONTEXT OF DRUG DISCOVERY**

Carla Newman⁽¹⁾ - **Andy West**⁽²⁾

⁽¹⁾ *The University of Nottingham, School of Pharmacy, Nottingham* - ⁽²⁾ *GSK, Bioimaging, Stevenage*

Keywords: LiveCell MS, intracellular drug, single cell metabolomics

Introduction

In the last few decades, the pharmaceutical industry has transformed people's lives. However, the development of new drugs is becoming increasingly and a paradigm shift in the drug discovery workflow is required to reduce attrition and transform conventional drug screening assays into translatable analytical techniques for the analysis of drugs in complex environments,

both in-vitro and ex-vivo. The ability to visualise unlabelled compounds inside the cell at physiological dosages can offer valuable insight into the compound behaviour both on and off-target.

Methods

LiveCell MS is a semi-automated methodology that allows the collection of intracellular content using a modified CQ1 imaging system developed by Yokowaga. The instrument is equipped with a confocal microscope that allowed bright field imaging as well as fluorescence imaging with 4 lasers (405, 488, 561 and 640 nm). Sampling was performed using the tips developed by Professor Masujima (reference below). The tip, holding the cellular contents, was then used for static nanospray of the contents into an Orbitrap Fusion Lumos (Thermo Scientific) and the resulting data processed using Compound Discoverer (Thermo Scientific).

Results

In this study, we show the applicability of the LiveCell MS technology to drug discovery, as it is crucial to identify compound and its metabolites when incubated in a mammalian cell at a therapeutic dose.

We report on the validation studies performed using the LiveCell MS platform. In the first of these studies we assess the dynamic range of the assay by directly comparing the results obtained with a state of the art metabolomics assay using bulk samples with our single cell technology.

In the second validation study we explore the effects of amiodarone, an anti-arrhythmia drug, in a macrophage model (rat-derived NR8383 cell line) identifying the drug intracellularly and the detectable endogenous metabolites.

Conclusions

We demonstrate that dosed compound can be identified in a single cell after sampling using the modified CQ1. Endogenous metabolites can also be identified, that can further understanding of the drug's mechanism. This technique has direct relevance for assessing compound effects on disease relevant cells and its low sample requirement makes it applicable to studying rare cell types. The use of high content imaging system enables the effect of compounds on live cells to be studied and suitable time points selected for sampling cell contents.

Novel Aspect

Identification of drug compounds and endogenous metabolites using an automated sampling technique and well as in combination with high-resolution mass spectrometry.

References

1. Fujii T., Matsuda S., Tejedor M.J., Esaki T., Sakane I., Mizuno H., Tsuyama N., Masujima T., Nature Protocols volume 10, pages 1445–1456 (2015) doi:10.1038/nprot.2015.084
2. Masujima T. Anal. Sci. 25, 953-60 (2009).
3. Mizuno H., Tsuyama N., Harada T., Masujima T., J. Mass Spectrom. 43, 1692-700 (2008) doi: 10.1002/jms.1460.
4. Guillaume-Gentil O., Rey O., Kiefer P., Ibáñez A.J., Steinhoff R., Brönnimann R., Dorwling-Carter L., Zambelli T., Zenobi R., Vorholt J.A., Anal. Chem. 89, 5017–5023 (2017) doi: 10.1021/acs.analchem.7b00367

TOr.02-4 IN SITU MASS SPECTROMETRY UNRAVELS METABOLIC DIFFERENCES BETWEEN CELLS OF THE LIVE VERTEBRATE EMBRYO (XENOPUS LAEVIS)

Erika Portero⁽¹⁾ - Sally Moody⁽²⁾ - Peter Nemes⁽¹⁾

⁽¹⁾ University of Maryland, Chemistry and Biochemistry, College Park - ⁽²⁾ The George Washington University, Anatomy and Regenerative Biology, Washington DC

Keywords: Single-cell, capillary electrophoresis, metabolomics, *Xenopus laevis*

Introduction

Metabolomics with single-cell resolution opens new frontiers in the study of cell heterogeneity [1] and establishment of functional differences during embryonic development [2]. Mass spectrometry (MS) is the technology of choice for metabolomics, but sensitivity limitations in detection hinder the application of traditional metabolomics to single embryonic cells. Here, we present a single-cell MS approach that enables the extraction and detection of metabolites from single cells from live *Xenopus laevis* embryos with high technical and biological reproducibility. We conduct *in situ*, minimally invasive sampling of identified embryonic cells to peer into metabolic pathways active during early embryogenesis. To enhance metabolite coverage, we integrate anionic and cationic analyses using a custom-built capillary electrophoresis electrospray ionization (CE-ESI) MS platform.

Methods

Single dorsal D1 and ventral V1 cells were identified in the 8-cell *Xenopus laevis* embryo using cell fate maps, pigmentation, and morphology criteria. Our CE-MS system [3] consisted of a microprobe to aspirate ~10 nL portion of the cell's cytoplasm, which was immediately extracted in 40 % acetonitrile/40% methanol and water. Metabolites were detected in a custom-built CE-ESI system coupled to a high-resolution quadrupole-time-of-flight mass spectrometer (Impact HD, Bruker).

Results

Our CE-ESI-MS platform provided 60 amol sensitivity and a quantitative range of 4–5 orders of magnitude. Trace-sensitive detection by this platform allowed us to detect ~250 cationic and ~150 anionic metabolite features in the positive and negative ion-mode. We identified 70 of the cationic features with high confidence as metabolites based on accurate mass, tandem MS, and migration time information against reference standards. Likewise, 50 of the anionic features were putatively assigned based on accurate mass match against data in Metlin [4]. These quantitative metadata enabled us to compare metabolic activities between cells along the dorsal-ventral developmental axis of the vertebrate embryo using multivariate and statistical analysis. Surprisingly, we observed clustering of metabolites based on their cell type of origin. Moreover, pathway analysis revealed the arginine-proline metabolism pathway (p-value 2.24×10^{-12}) and the glycine-serine-threonine metabolism pathway (p-value 5.13×10^{-10}) to be of significance, suggesting the importance of these pathways during early development.

Conclusions

We developed a unique single-cell MS technology that enabled the metabolic characterization of identified single cells in the live vertebrate embryo of *Xenopus laevis*. Quantitative data from this study revealed previously unknown metabolic cell heterogeneity between the dorsal-ventral axis of the 8-cell embryo, raising the possibility that metabolites also contribute to cell fate specification.

Novel Aspect

Enabling *in situ* sampling and detection of cationic and anionic metabolites from single-cells during early development to uncover cell heterogeneity.

References

1. Zenobi, R., *Science*, 342,1243259 (2013).
2. Fessenden, M., *Nature*, 540, 153-155 (2016).

3. Onjiko, R.M., Portero, E.P., Moody, S.A., Nemes, P., *Anal. Chem.*, 89, 7069 (2017).
4. Metlin: The Scripps Research Institute, <https://metlin.scripps.edu/>

TOr.03 - Noncovalent interactions

Chair: Renato Zenobi

TOr.03-1 **Keynote: MEMBRANE PROTEINS - THE LIPID CONNECTION**

Carol V. Robinson

Physical and Theoretical Chemistry Laboratory, University of Oxford

Keywords: Membrane proteins, detergent micelles, G protein coupled receptors

Introduction

Mass spectrometry (MS) has the potential to monitor different modes of lipid binding to membrane protein complexes. Initial studies monitored the addition of lipids and deduced the kinetic and thermodynamic effects of lipid binding to proteins. Recently we have identified lipids already present, explicitly in plugs, annular rings or cavities. These lipids have higher residence times and consequently can be quantified and characterized by mass spectrometry.

Methods

We have used native MS of membrane proteins liberated from micelles to study protein lipid interactions. Using an Orbitrap MS platform developed for membrane proteins (1) we show that it is possible to distinguish lipids, cofactors and drugs while in direct contact with the membrane protein target.

Results

While most steps of lipid II transfer across the cytoplasmic membrane are understood, identifying the lipid II flippase has yielded conflicting results. We used native mass spectrometry to characterise lipid II binding to two candidate proteins MurJ and FtsW. We showed that lipid II binding competes with cardiolipin to reduce lipid II binding to MurJ (2). We also showed how drug binding mediates phosphorylation in a class A GPCR (P2Y₁R) (3) and that purification of three further class A GPCRs results in binding of specific lipids. These lipids in turn affect downstream coupling to certain types of G-proteins and consequently appear to enhance selective coupling primarily to class A GPCRs.

Conclusions

The ability to detect lipid binding which have both structural and functional effects, for example on lipid flippase activity or G-protein coupling, is a powerful and important attribute of the study of membrane proteins by mass spectrometry.

Novel Aspect

First mass spectra of folded GPCR complexes and demonstration of competition between lipids affecting function.

References

1. Gault, J., Donlan, J. A., Liko, I., Hopper, J. T., Gupta, K., Housden, N. G., Struwe, W. B., Marty, M. T., Mize, T., Bechara, C., Zhu, Y., Wu, B., Kleanthous, C., Belov, M., Damoc, E., Makarov, A., and Robinson, C. V. *Nature Methods* 13, 333-336 (2016) .

- Bolla, J. R., Sauer, J. B., Wu, D., Mehmood, S., Allison, T. M., and Robinson, C. V. Direct observation of the influence of cardiolipin and antibiotics on lipid II binding to MurJ. *Nature Chemistry* 10, 363-371. (2018).
- Yen, H. Y., Hopper, J. T. S., Liko, I., Allison, T. M., Zhu, Y., Wang, D., Stegmann, M., Mohammed, S., Wu, B., and Robinson, C. V. (2017) Ligand binding to a G protein-coupled receptor captured in a mass spectrometer. *Science Advances* 3, e1701016 (2017).

TOr.03-2 ARYL BIS-SULFONAMIDES BIND TO THE ACTIVE SITE OF A HOMOTRIMERIC ISOPRENOID BIOSYNTHESIS ENZYME ISPF AND EXTRACT THE ESSENTIAL DIVALENT METAL CATION COFACTOR

Katharina Root⁽¹⁾ - **Konstantin Barylyuk**⁽²⁾ - **Anatol Schwab**⁽¹⁾ - **Jonas Thelemann**⁽¹⁾ - **Boris Illarionov**⁽³⁾ - **Markus Fischer**⁽³⁾ - **François Diederich**⁽⁴⁾ - **Renato Zenobi**⁽¹⁾

ETH Zurich, Department of Chemistry and Applied Biosciences, Zurich⁽¹⁾ - ⁽²⁾ University of Cambridge, Department of Biochemistry, Cambridge - ⁽³⁾ University of Hamburg, School of Food Science, Hamburg

Keywords: Native ESI-MS, oligomeric enzyme, binding mechanism, aryl bis-sulfonamides, IspF

Introduction

Characterizing the mode of action of non-covalent inhibitors in multisubunit enzymes often presents a great challenge. Most of the conventional used techniques often fail because they do not reveal how the distribution of individual ligand-bound states changes with ligand concentration but instead describe an average property of the ensemble as a function of ligand concentration. Thus, mixed models comprising both competitive and noncompetitive binding are often used that proved very little insight into the mechanisms of ligand binding and recognition.

Methods

Native electrospray ionization mass spectrometry (ESI-MS) has recently emerged as a powerful technique for direct, label-free analysis of intact protein complexes. Native ESI-MS has been increasingly recognized for its ability to directly probe such properties of intact macromolecular assemblies as the protein-ligand complex stoichiometry, binding affinities, and allosteric effects.

Results

Here we apply native ESI-MS to comprehensively characterize an isoprenoid biosynthesis enzyme IspF from *Arabidopsis thaliana*, an example of homomeric protein complex with multiple binding sites for several types of ligands, including a metal cofactor and a synthetic inhibitor. Based on these combined with computational modeling, we propose a mechanism of AtIspF inhibition by aryl bis-sulfonamides that involves both the competition with the substrate for the ligand-binding pocket and the extraction of Zn²⁺ from the enzyme active site. This inhibition mode is therefore mixed competitive and non-competitive, the latter exerting a key inhibitory effect on the enzyme activity.

Conclusions

The results of our study deliver a profound insight into the mechanisms of AtIspF action and inhibition, open new perspectives for designing inhibitors of this important drug target, and demonstrate the applicability and value of the native ESI-MS approach for deep analysis of complex biomolecular binding equilibria.

Novel Aspect

While standard biophysical techniques failed to reveal the mode of action of aryl-sulfonamide-based inhibitors of AtIspF, direct native ESI-MS titrations of the protein with the ligands and ligand competition assays allowed us to accurately capture the solution-phase protein-ligand binding equilibria in full complexity and detail.

TOr.03-3 NATIVE MS OF NONCOVALENT RNA:PROTEIN COMPLEXES

Vicki Wysocki

Dept. of Chemistry and Biochemistry, Ohio State University

Keywords: Native MS; protein complex; nucleoprotein; surface-induced dissociation; structural biology

Introduction

Characterization of the overall topology and inter-subunit contacts of non-covalent protein complexes, and their assembly/disassembly and unfolding pathways, is critical because protein complexes regulate key biological processes, including those important in understanding and controlling disease. Native mass spectrometry is an approach that provides critical structural information at high throughput on low sample amounts.

Methods

The power of native MS increases when coupled to ion mobility (IM-MS), a technique that measures rotationally averaged collisional cross sections and thus information on conformational changes, or to high resolution mass spectrometry (HRMS). This presentation illustrates CID, SID/IM, and SID HRMS for characterization of topology, inter subunit connectivity, and other structural features of multimeric nucleoprotein complexes.

Results

Native MS is being used in our lab for characterization of soluble, membrane, and nucleoprotein complexes. Native MS results for multiple systems will be presented in this talk, with a primary focus on RNA:protein complexes. The first nucleoprotein system under investigation (with Venkat Gopalan of OSU) is RNase P, an endonuclease that catalyzes 5' maturation of t-RNAs. While we've been able to bind four RNase P proteins to the RNA or a fifth protein to kink turns in the RNA, we are still optimizing conditions to simultaneously bind all five proteins. A second system under investigation (with Karin Musier-Forsyth of OSU) involves the retroviral Gag protein involved in HIV-1 genome packaging and viral assembly, a system for which we've measured 3:1 Gag:Psi complexes. A third nucleoprotein project (with Sarah Woodson of Johns Hopkins) involves the hexameric RNA chaperone Hfq bound to RNA and recruiting messenger RNA, a system for which surface-induced dissociation of the hexameric ring is altered by the presence and absence of the bound RNA.

Conclusions

Native MS coupled to HCD, SID, ion mobility, or high resolution mass spectrometry, along with complementary tools such as covalent labeling MS, provides information on conditions needed for optimal formation of the protein:RNA complexes (ionic strength, Mg²⁺ concentration, protein:RNA ratios) and on the protein-RNA connectivity and structure. Experiments are underway to determine the best ways to calibrate the CCSs of the RNAs and RNA:protein complexes.

Novel Aspect

Native MS and SID of nucleoprotein complexes provides structural information unavailable or difficult to obtain by other structural biology tools.

TOr.03-4 EXPLORING THE INHIBITION MECHANISMS OF AMYLOID PROTOFIBRILS USING FOURIER TRANSFORM ION CYCLOTRON RESONANCE-MASS SPECTROMETRY (FTICR-MS)

Pui Yiu Yuko Lam⁽¹⁾ - Cookson K. C. Chiu⁽¹⁾ - Christopher A. Wootton⁽¹⁾ - Ji-Inn Song⁽¹⁾ - Meng Li⁽¹⁾ - Ian Hands-Portman⁽²⁾ - Mark P. Barrow⁽¹⁾ - Peter B. O'Connor⁽¹⁾

(¹) University of Warwick, Department of Chemistry, Coventry - (²) University of Warwick, Department of Life Science, Coventry

Keywords: Amyloid protein, Inhibition mechanism, Deamidation, Non-covalent interaction, Aggregation

Introduction

Human Islet Amyloid Polypeptide (hIAPP) is an amyloidogenic protein that aggregates rapidly in humans.^[1-5] Classic structure-based design for therapeutics cannot apply to amyloid protein drug development as most amyloid proteins are inherently disordered.^[6-7] Herein, Fourier Transform Ion Cyclotron Resonance (FTICR) MS was used to study the potential inhibitors, identify binding regions through top-down MS/MS, and correlate with fluorescence and transmission electron microscopy (TEM) measurements to understand target regions for aggregation inhibition mechanisms.

Methods

Collisionally activated dissociation (CAD) MS/MS and electron capture dissociation (ECD) MS/MS in a 12T FTICR MS were used to show the interaction strengths and regions between hIAPP/deamidated hIAPP and the potential inhibitors respectively.

The performance of the aggregation inhibitors on preventing hIAPP/ deamidated hIAPP aggregation was quantified by a fluorescence spectrometer and relative MS quantification to determine the amount of fibrils formed in the incubated solutions. The overall fibrils structures were observed by TEM.

Results

Insulin, 3-APS, BISA, EGCG are all reported inhibitors of hIAPP aggregation. Insulin was found to specifically bind to the critical aggregation region of hIAPP (Ser-29 to Asn-35) which inhibited the formation of the early oligomers and the mature fibrils. The hIAPP in the heterodimer complex was easily converted back to soluble hIAPP in an aqueous solution, thus, 100% recovery of the soluble hIAPP was observed in an incubated solution.

On the contrary, EGCG bound non-specifically to hIAPP and the interaction pattern of EGCG was very similar to the potential drugs which were found to be ineffective at inhibiting hIAPP aggregation (3-APS and BISA). Instead, EGCG was found to aggregate rapidly with hIAPP in a non-fibrillary complex to inhibit the aggregation.

For deamidated hIAPP, CAD MS/MS results showed insulin had a weaker interaction with the deamidated hIAPP compared to the wild-type and mature amorphous aggregates were formed during the experiment.

Conclusions

The results demonstrated there are two potential inhibition mechanisms adopted by the potential inhibitors to prevent hIAPP aggregation. Insulin directly inhibits hIAPP fibrillization by specific interaction with target sequence regions, and EGCG rapidly redirects hIAPP into non-toxic amorphous aggregates which prevents the amyloid fibrils formation. Only EGCG inhibitors with redirecting the protein into non-toxic amorphous aggregates was found to be applicable in inhibiting the deamidated hIAPP aggregation.

Novel Aspect

Understanding underlying mechanisms that prevent amyloid aggregation are the key to amyloid drug development

References

1. Westermark, P.; Wernstedt, C.; Wilander, E.; Hayden, D. W.; O'Brien, T. D.; Johnson, K. H., *Proc. Natl. Acad. Sci. USA* 84 (11), 3881-3885 (1987).

- Cooper, G.; Willis, A.; Clark, A.; Turner, R.; Sim, R.; Reid, K., *Proc. Natl. Acad. Sci. USA* 84 (23), 8628-8632 (1987).
- Luca, S.; Yau, W. M.; Leapman, R.; Tycko, R., *J. Biochem.* 46 (47), 13505-13522 (2007).
- Cao, P.; Marek, P.; Noor, H.; Patsalo, V.; Tu, L.-H.; Wang, H.; Abedini, A.; Raleigh, D. P., *FEBS Lett.* 587 (8), 1106-1118 (2013).
- Akter, R.; Cao, P.; Noor, H.; Ridgway, Z.; Tu, L. H.; Wang, H.; Wong, A. G.; Zhang, X. X.; Abedini, A.; Schmidt, A. M.; Raleigh, D. P., *J. Diabetes Res.*, (2016).
- Noble, M. E.; Endicott, J. A.; Johnson, L. N., *Science* 303 (5665), 1800-1805 (2004).
- Hajduk, P. J.; Greer, J., *Nat. Rev. Drug Discovery* 6 (3), 211-219 (2007).

TO.r.03-5 **DIFFERENT MEANS OF SOLUBILIZATION OF MEMBRANE PROTEINS FOR NATIVE MS**

Nina Morgner - Oliver Peetz - Nils Hellwig

Frankfurt University, Frankfurt University / Physical Chemistry, Frankfurt

Keywords: membrane proteins, nanodiscs, SMALPs, native MS, LILBID

Introduction

Membrane proteins are of high interest, but still underrepresented in research due to the inherent difficulties for many investigation methods, arising from their hydrophobic nature. Depending on the feature of the membrane protein complex the researcher is interested in (protein stoichiometry, lipid affinity, annular lipids...), different means of solubilisation can enable MS analysis.

Methods

LILBID (Laser Induced Liquid Bead Ion Desorption) is an ionization method which employs a droplet generator to produce analyte droplets of 30-50 µm diameter at a frequency of 10 Hz [1]. These droplets are irradiated by a mid-IR laser leading to the explosive expansion of the droplet. The solvated ions are set free and are mass analyzed.

Results

We investigated the usability of artificial membranes for the study of membrane protein complexes with LILBID-MS. Lipid bilayer mimics as nanodiscs and liposomes allow the analysis of the oligomerization state of a protein complex. Cell-free expressed membrane proteins can be solvated by lipid bilayer mimics without ever being in contact with detergents. This allows for example the study of lipid dependent oligomerization [2]. Even SMALPs, which allow to cut protein complexes out of living cells preserving their annular lipid belt, can be used to determine the constituting proteins of a complex, in dependence of growth conditions of the cells.

Conclusions

LILBID-MS is well suited for the analysis of membrane proteins complexes solubilized by detergent or different lipid bilayer mimics, such as nanodiscs, liposomes or SMALPs. This allows analysis of the complexes' constituting proteins, as well as their dependence on specific lipids or ligands.

Novel Aspect

Analysis of membrane protein complexes directly out of lipid bilayer mimics

References

- Morgner N., Barth H. D., Brutschy B., *Australian Journal of Chemistry*, 59 (2), 109-114 (2006).
- Henrich E., Peetz O., Hein C., Laguerre A., Hoffmann, Hoffmann J., Dötsch V., Bernhard F., Morgner N., *eLife*, 6, e20954 (2017)

TOr.04 - High resolution MS: instrumentation, methods & applications

Chairs: Evgeny Nikolaev, Gianluca Giorgi

TOr.04-1 **Keynote: LIFE AT THE FRONTIER: 21 TESLA FOURIER TRANSFORM ION CYCLOTRON RESONANCE MASS SPECTROMETRY**

Alan G. Marshall⁽¹⁾ - **Martha L. Chacón-Patiño**⁽²⁾ - **Yuri E. Corilo**⁽²⁾ - **Christopher L. Hendrickson**⁽²⁾ - **Lidong He**⁽¹⁾ - **Sydney F. Niles**⁽¹⁾ - **Jonathan C. Putman**⁽¹⁾ - **Steven M. Rowland**⁽³⁾ - **Donald F. Smith**⁽²⁾ - **Rebecca L. Ware**⁽¹⁾ - **Chad R. Weisbrod**⁽²⁾ - **Ryan P. Rodgers**⁽²⁾

⁽¹⁾ Department of Chemistry & Biochemistry, Florida State University, 95 Chieftain Way, Tallahassee, FL 32306 - ⁽²⁾ Ion Cyclotron Resonance Program, National High Magnetic Field Laboratory, 1800 East Paul Dirac Drive, Tallahassee, FL 32310 - ⁽³⁾ Florida State University Future Fuels Institute, 1800 E. Paul Dirac Dr., Tallahassee, FL 32310

Keywords: 21 tesla FT-ICR MS, top-down MS/MS, monoclonal antibodies, hemoglobin variants, petroleomics

Introduction

FT-ICR MS offers at least 10X times higher mass resolving power than any other mass analyzer, and is thus the mass analyzer of choice for “omics” applications. 21 tesla FT-ICR MS [1] enables ultrahigh mass resolving power at on-line HPLC speed (< 1 second per acquired transient) for high-throughput analysis of complex organic mixtures: e.g. assignment of >450,000 distinct elemental compositions from a single LC FT-ICR mass chromatogram of Atjabasca bitumen.

Methods

We perform protein top-down and middle-down MS/MS (front-end ETD [2], CID, or UVPD; in-cell IRMPD). Our ICR cell is dynamically harmonized, with 120° azimuthal excitation/detection (for improved S/N and elimination of 3rd harmonics) and segmented electrodes (for 3X ICR frequency detection, and thus mass resolving power equivalent to 3X magnetic field strength). Broadband phase correction yields absorption-mode mass spectra (for 2X LC/MS mass resolving power).

Results

Based on various front-end extraction/chromatography procedures, we identify up to 10X more elemental compositions in petroleum than from any single broadband mass spectrum. For intact proteins, we report monoisotopic mass measurements 10x more accurate than commercial averaging-based methods. For example, we identify the class, sequence variants, and multiple glycoforms of monoclonal antibodies from multiple myeloma patient serum, with 100% match to DNA germline sequencing, and without false negatives due to incomplete bone marrow aspiration. For human hemoglobin, we are able to identify sequence variants as well as thalassemia (from, e.g., the abundance ratio of delta chain/beta chain). Work supported by NSF Division of Materials Research through NSF DMR-11-57490 and DMR-1644779, The Gulf of Mexico Research Initiative to the Deep-C Consortium, the Florida State University Future Fuels Institute, and the State of Florida.

Conclusions

FT-ICR mass resolving power and data acquisition speed increase directly with magnetic field strength, B_0 . Also, mass accuracy, dynamic range, and peak non-coalescence scale as B_0^2 . This talk will present new benchmarks for FT-ICR MS performance, based not only on 21 tesla magnetic field (highest in the world), but also on recent advances in front-end separations and MS/MS, ICR excitation and detection, and data reduction (mass calibration; phase correction).

Novel Aspect

21 tesla FT-ICR MS [1], with front-end ETD [2], dynamically harmonized ICR cell with 120° azimuthal excitation/detection and tripled frequency detection

References

1. Hendrickson, C. L.; Quinn, J. P.; Kaiser, N. K.; Smith, D. F.; Blakney, G. T.; Chen, T.; Marshall, A. G.; Weisbrod, C. R.; Beu, S. C. *J. Am. Soc. Mass Spectrom.*, 26, 1626-1632 (2015)
2. Weisbrod, C. R.; Kaiser, N. K.; Syka, J. E. P.; Early, D.; Mullen, C.; Dunyach, J.-J.; English, A. M.; Anderson, L. C.; Blakney, G. T.; Shabanowitz, J.; Hendrickson C. L.; Marshall, A. G.; Hunt, D. F. *J. Am. Soc. Mass Spectrom.*, 28, 1787-1795 (2017)

TO.r.04-2 BRIDGING THE GAP BETWEEN IMAGES AND MOLECULES: PARALLEL ACQUISITION OF FTMS MASS SPECTROMETRY IMAGING AND MS/MS DATA COUPLED WITH AUTOMATED AND HIGH CONFIDENCE LIPID IDENTIFICATIONS

Shane Ellis⁽¹⁾ - Martin Paine⁽¹⁾ - Gert Eijkel⁽¹⁾ - Josch Pauling⁽²⁾ - Peter Husen⁽³⁾ - Mark Jervelund⁽²⁾ - Martin Hermansson⁽²⁾ - Christer Ejsing⁽²⁾ - Ron M. A. Heeren⁽¹⁾

⁽¹⁾Maastricht University, The Maastricht Multimodal Molecular Imaging Institute (M4I), Maastricht - ⁽²⁾University of Southern Denmark, Department of Biochemistry and Molecular Biology, Villum Center for Bioanalytical Sciences, Odense - ⁽³⁾University of Southern Denmark, Department of Physics, Chemistry and Pharmacy, Odense

Keywords: mass spectrometry imaging; lipids; tandem mass spectrometry; Fourier transform mass spectrometry; automated identification

Introduction

Mass spectrometry imaging (MSI) enables visualizations of lipid distributions throughout tissues. Yet, it has struggled to translate images to underlying lipid biochemistry. A significant contributor is the lack of automated and broad MS/MS acquisition and structural identifications tool. We address this limitation via an innovative MSI approach enabling parallel acquisition of MS/MS and FTMS-MSI of every ion without prolonging MSI acquisition time.

Methods

Data were acquired using an Orbitrap Elite and a dual MALDI/ESI source (Spectrograph, LLC). Norharmane matrix and a 20x40 µm step size were used for MSI [1]. For DDA-Imaging IT-MS/MS were acquired in parallel with each Orbitrap transient. Acquisition of MS/MS for every *m/z* was achieved without prolonging the FTMS-MSI acquisition. ALEX¹²³ software was used to assign both sum-composition and molecular-lipid identifications utilising both FTMS and IT-MS/MS data [2].

Results

Exploiting the parallel operation of the Orbitrap Elite our DDA-Imaging-method yielded both high mass resolution FTMS (encoding both the MSI data and sum-composition level lipid assignments) and high quality IT-MS/MS of nearly every ion observed in the MSI experiment. Both a 240,000 resolution FTMS scan and an IT-MS/MS spectrum could be obtained in parallel from each 40 µm region of tissue. Even for low abundance lipids having relative intensities less than 1% high contrast images and high quality MS/MS data were acquired.

We used ALEX¹²³ to perform automated identification of detected lipids using both the high mass accuracy (+/- 2 ppm tolerance) and IT-MS/MS data. We automatically identified 165 unique sum-composition-level lipid species and 113 unique molecular lipid species that were present in all replicates of rat cerebellum tissue. All IDs are supported by both FTMS and MS/MS and in relevant cases, cross-correlation of positive- and negative-ion mode data. For all lipid species their spatial distribution is also acquired via the FTMS data.

Conclusions

Our method overcomes the current limitations of broadband molecular identifications in MSI. By exploiting a hybrid ion trap/Orbitrap instrument, MS/MS data area acquired without expanding additional experiment. To fully exploit the information-rich datasets, ALEX¹²³ software provides automated and high confidence lipid identifications in parallel to MSI experiments, allowing broad assignment of MS-images to identified lipids.

Novel Aspect

Parallel FTMS-MSI and MS/MS combined with ALEX¹²³ software provides both MSI and broad structural identification data in a single acquisition.

References

1. Belov, M.E. et al. *Anal. Chem.* **89**, 7493-7501 (2017).
2. Almeida, R. et al. *J. Am. Soc. Mass. Spectrom.* **26**, 133-148 (2015).

TOr.04-3 **LARGE-SCALE RING TRIAL METABOLOMICS STUDY OF HUMAN BLOOD PLASMA USING DIVERSE FT-ICR-MS SYSTEMS**

Sara Forcisi^(1,2) - **Franco Moritz**⁽¹⁾ - **Basem Kanawati**⁽¹⁾ - **Jenny Uhl**⁽¹⁾ - **Christopher Thompson**⁽³⁾ - **Aiko Barsch**⁽⁴⁾ - **Gary Kruppa**⁽³⁾ - **Michael Easterling**⁽³⁾ - **Ring Trial Consortium**⁽¹⁾ - **Philippe Schmitt-Kopplin**^(1,2,5)

⁽¹⁾ Research Unit Analytical BioGeoChemistry, Department of Environmental Sciences, Helmholtz Zentrum München, Ingolstädter Landstraße 1, 85764 Neuherberg, Germany - ⁽²⁾ German Center for Diabetes Research (DZD), Ingolstädter Landstraße 1, 85764 Neuherberg, Germany - ⁽³⁾ Bruker Daltonics Inc., 40 Manning Road, Billerica, MA 01821, USA - ⁽⁴⁾ Bruker Daltonics GmbH, Life Sciences, Fahrenheitstrasse 4, 28359 Bremen, Germany - ⁽⁵⁾ Chair of Analytical Food Chemistry, Technische Universität München, Alte Akademie 10, 85354 Freising-Weihenstephan, Germany

Keywords: FT-ICR-MS, Ring Trial, human blood plasma, metabolomics

Introduction

A ring trial study of main blood plasma metabolites, in collaboration with over 17 different FT-ICR-MS laboratories worldwide, was organized and administrated at the Helmholtz Center Munich in a strong cooperation with Bruker Daltonics. The aim of this study is to investigate and to validate the implementation of several important ultra-high resolution FT-ICR-MS instruments for performing phenomics studies.

Methods

A novel standard operating procedure (SOP) was developed at the Helmholtz Center Munich, to assess deep metabolic information on human blood plasma samples[1]. Spiking studies with mixtures of chemical standards were performed to challenge the instruments' performance.

Results

The overlap of known and unknown metabolic compounds (features), reproducibly detected in all labs, is highlighted. In addition, key instrumental specificities between different FT-ICR-MS mass spectrometers are investigated, such as reproducibility, resolution, FWHM, peak shape, measurement sensitivity, dynamic range and mass measurement accuracy.

Conclusions

The study explores the capabilities of different FT-ICR-MS instruments with diverse ICR cells working under different magnetic field strengths, and with different operational cyclotron frequencies (1st versus 2nd ion frequency measurements).

Novel Aspect

This study provides an unprecedented large-scale investigation to reliably perform accurate mass measurements for serving phenomics studies in the future.

References

1. Forcisi S., Moritz F., Lucio M., Lehmann R., Stefan N., Schmitt-Kopplin P., *Anal. Chem.*, 87, 8917-24 (2015).

TO.r.04-4 SINGLE CELL METABOLOMICS WITH LAESI-21T-FTICR-MS

Ljiljana Paša-Tolić⁽¹⁾ - Christopher R. Anderton⁽¹⁾ - Dusan Veličković⁽¹⁾ - Sylwia A. Stopka⁽²⁾ - Laith Samarah⁽²⁾ - Rikkita Khattar⁽²⁾ - Jared B. Shaw⁽¹⁾ - Beverly J. Agtuca⁽³⁾ - Gary Stacey⁽³⁾ - Akos Vertes⁽²⁾

⁽¹⁾ Pacific Northwest National Lab, Environmental Molecular Sciences Lab, Richland - ⁽²⁾ The George Washington University, Department of Chemistry, Washington - ⁽³⁾ University of Missouri, C. S. Bond Life Sciences Center, Columbia

Keywords: MSI, FTICR, LAESI, plant-microbe

Introduction

The ability to visualize complex metabolic behavior of living organisms *in situ* and at the single cell level represents a bioanalytical grand challenge that has yet to be fully realized. Herein, we present the combination of fiber-LAESI with fluorescence microscopy and ultrahigh resolution 21T-FTICR-MS to address these needs.

Methods

To provide simultaneous spatial distributions for hundreds of metabolites in biological tissues, an LAESI ion source was coupled with a 21T FTICR-MS system. Utilizing the ultra-high mass resolving power, isotopic fine structure readout was used for metabolite annotations. Single-cell spatial resolution was achieved by delivering mid-IR laser pulses through a sharpened optical fiber, while a dual modality microscopy allowed for precise targeting of cells.

Results

Spatial heterogeneity was tackled using LAESI coupled with the 21T-FTICR-MS. This unique instrumental configuration provides exquisite mass and spatial resolution, down to the single cell level. A dual channel microscope capable of both simultaneous brightfield and fluorescence imaging was combined with the fiber-LAESI-MS to enable selection of specific cells. For instance, cells in a soybean root nodule infected by GFP labeled soil bacteria were targeted for analysis in a heterogeneous population of infected and uninfected cells. The superior mass resolution and accuracy offered by the 21T FTICR-MS facilitated identification of differentially abundant metabolites in soybean root nodule cells infected by *wt* and *nifH*- rhizobia. Notably, this novel capability facilitated characterization of organometallic molecules due to their natural isotope patterns (e.g. ⁵⁷Fe isotopologues). Preliminary data indicates significantly lower abundance of heme b, disaccharide, and S-adenosylmethionine in the *nifH*- mutant.

Conclusions

We demonstrate that LAESI coupled with the 21T FTICR-MS holds tremendous potential for *in-situ* single cell metabolomics as illustrated by the study of plant-microbe interactions.

Novel Aspect

In situ single cell analysis by fiber-LAESI-MS promises new insights into cellular heterogeneity and metabolic noise.

TOr.04-5 **COMPREHENSIVE HIGH-RESOLUTION MASS SPECTROMETRIC EVOLVED GAS ANALYSIS IN THE CONTEXT OF PETROLEOMICS: ANALYSIS OF BITUMEN, HEAVY CRUDE OILS AND ASPHALTENES**

Christopher Ruger⁽¹⁾ - **Anika Neumann**⁽¹⁾ - **Christoph Grimmer**⁽¹⁾ - **Uwe Kafer**⁽²⁾ - **Martin Sklorz**⁽¹⁾ - **Thomas Groger**⁽²⁾ - **Thorsten Streibel**⁽¹⁾ - **Ralf Zimmermann**⁽¹⁾

⁽¹⁾University of Rostock, Institute of Chemistry, Rostock - ⁽²⁾Helmholtz Zentrum Munchen, Comprehensive Molecular Analytics, Munchen

Keywords: evolved gas analysis, high resolution mass spectrometry, Petroleomics, atmospheric pressure ionization, photo ionization

Introduction

Petroleomic materials remain to be one of the greatest challenges in chemical analysis. The conversion of heavy oils or exploration of alternative sources raises new challenging issues, such as the formation of deposits or efficient processing of residues. Thermal analysis using evolved gas analysis coupled to mass spectrometry is one powerful approach for deciphering these materials.

Methods

For this purpose, several thermal analysis systems coupled to different mass spectrometric analyzers, such as ultra-high resolution FT-MS [1,2] and TOF-MS, were applied to study low-viscous and solid petroleomic materials. The setups varied in ionization techniques covering soft vacuum and atmospheric pressure approaches, and hard electron impact allowing in combination an in depth-chemical description of the desorbed and pyrolyzed material. [3]

Results

Besides the description of heavy oils and bitumen, asphaltenes were in the primary focus of the study. Generally, for the intended pyrolysis of the asphaltenes, an ultra-complex molecular pattern is revealed. This molecular space could be comprehensively analyzed by ultra-high resolution mass spectrometry equipped with atmospheric pressure chemical and photo ionization revealing the highly aromatic and Sulphur-containing core structures, whereas vacuum single photon ionization was able to analyze the thermally fragmented side chain distribution. Moreover, a vacuum direct inlet probe equipped with electron ionization was able to broaden the picture by evaporating larger components prior thermal decomposition and thereby allow for additional structural information.

Conclusions

Summarizing, mass spectrometric evolved gas analysis has shown to be an efficient technique for Petroleomics. Most importantly the data combination of various setups enabled an in-depth chemical and structural analysis and hypothesizing of macromolecular structural features. Future work will focus on data integration as well as on structural information by adding a chromatographic separation between thermal release and mass spectrometric detection.

Novel Aspect

Data integration of several mass spectrometric evolved gas analyzer allowing comprehensive chemical analysis and structural motives of complex petrochemical materials.

References

1. Ruger, Miersch, Schwemer, Sklorz, Zimmermann, *Analytical Chemistry*, 87, 6493–6499 (2015).
2. Ruger, Neumann, Sklorz, Schwemer, Zimmermann, *Energy and Fuels*, 31, 13144–13158 (2017).

3. Ruger, Grimmer, Sklorz, Neumann, Streibel, Zimmermann, *Energy and Fuels*, 32, 2699–2711 (2018).

TOr.05 - MS imaging: instrumentation

Chairs: Liam McDonnell, Fulvio Magni

TOr.05-1 **Keynote: MOLECULAR IMAGING AT THE INTERFACE OF HEALTH AND DISEASE**

Ron M. A. Heeren

M4I, The Maastricht MultiModal Molecular Imaging institute, Maastricht University, Universiteitssingel, 50, 6229 ER Maastricht, The Netherlands

Keywords: Imaging MS, MALDI, SIMS, Clinical diagnostics, Isomers

Introduction

A comprehensive understanding of molecular patterns of health and disease is needed to pave the way for personalized medicine and tissue regeneration [1]. The best way to capture disease complexity is to chart and connect multilevel molecular information within a tissue, using mass spectrometry. Charting this territory with molecular maps from cells and tissue has become reality through the clinical implementation of imaging mass spectrometry [2]

Methods

High resolution FTMS technologies (FTICR-MS and Orbitrap MS) combined with smart ion chemistry, stable isotope labeling approaches, ambient ionization, new data acquisition approaches and funnel based MALDI ion sources allows us to address some of the open challenges that still exist in the field of imaging MS. In concert with high throughput MALDI-TOF imaging technologies truly translation clinical studies have become reality.

Results

We pushed the limit of spatial resolution, molecular resolution and maybe more importantly, structural resolution forward in imaging MS. A novel elevated pressure MALDI imaging ion source is employed to reveal local biochemistry at high throughput. In addition, we have employed OzID in combination with MALDI-MSI to reveal the distinctive regulatory role of lipid isomeric forms in a variety of diseases. These methods are employed for the generation of molecular atlases from cells and tissue and complemented with high throughput “omics” approaches. Multimodal approaches enable the study of clinical samples at a variety of molecular and spatial scales. The combination with tools from structural biology makes it possible to perform mass spectrometric experiments at length scales from cells to patients.

Conclusions

We have demonstrated how new MS based chemical microscopes target biomedical tissue analysis in various diseases as well as other chemically complex surfaces. In concert they elucidate the way in which local environments can influence molecular signaling pathways on various scales. The integration of this pathway information in a surgical setting is imminent, but innovations that push the boundaries of the technology and its application are still needed.

Novel Aspect

Imaging MS reveals new structural detail in tissue analysis.

References

1. Vaysse P.M., Heeren R.M.A., Porta T., Balluff B. *Analyst*, 142, 2690-2712 (2017)
2. Chughtai K., Heeren R.M.A. *Chem. Rev.*, 110, 3237-3277 (2010)

TO.r.05-2 **NOVEL HYPERSPECTRAL IMAGING MASS SPECTROMETRY TO DELINEATE AMYLOID POLYMORPHISM**

Jörg Hanrieder

University of Gothenburg, Department of Neuroscience and Physiology, Gothenburg

Keywords: MALDI imaging mass spectrometry, hyperspectral chemical imaging, electrooptic fluorescent probes, Alzheimer's disease, plaque polymorphism

Introduction

Amyloidogenic aggregation of beta-amyloid peptides into senile deposits is the major pathological hallmark of Alzheimer's disease (AD). While the exact mechanisms of AD pathogenesis are not fully understood, plaque pathology has been identified as critical, driving AD pathogenesis. Delineating plaque pathology in situ is however challenging and calls for novel, multimodal chemical imaging tools such as imaging mass spectrometry¹ in combination with other biochemical imaging techniques^{2,3}.

Methods

A hyperspectral chemical imaging paradigm is presented, implementing multimodal MALDI imaging MS of peptides and lipids³ along with chemical imaging using electrooptic fluorescent probes that recognize structural characteristics of A β aggregation dynamics⁴. The method was employed for probing A β plaque polymorphism in human AD brain as well as in transgenic AD mouse models⁵.

Results

The results allowed to identify various A β and lipid species associated with AD pathogenesis. Here, the data on in human and mouse tissue show both a significant increase of A β 1-40 in senile, neurotoxic plaques as compared to plaques of diffuse, pre-amyloid morphology, where no A β 1-40 was observed both in human AD as well as in control patients that exhibit amyloid pathology but no cognitive decline. Additional data in mice show that this increase in A β 1-40 was associated with plaque maturation over time, where both senile plaques and cerebrovascular deposits at 18 months contained higher A β 1-40/1-42 levels than at 12 months. A proportionally higher A β 1-42 content appears therefore to be characteristic for pre-mature deposits suggesting diffuse plaques to be precursors of senile plaques and that maturation into senile plaques is associated with AD pathogenesis. These data show for the first time what A β and lipid profiles are associated with differences in plaque morphology in progressing A β pathology.

Conclusions

Taken together, the here presented hyperspectral chemical imaging paradigm of individual, chemically and structurally delineated A β deposits provides insight in the chemical aspects underlying plaque polymorphism observed in Alzheimer's disease pathology that is not facilitated by other biochemical techniques. The results show that it is possible to attribute structural features of plaques to A β profiles and highlight the role of different A β isoforms and associated neuronal lipids in formation of distinct biochemical conformation of A β aggregates.

Novel Aspect

The data highlight the potential of implanting different MALDI imaging modes with other chemical imaging tools for comprehensively probing molecular disease pathology.

References

1. McDonnell, L. A.; Heeren, R. M. A. *Mass Spectrometry Reviews* 26, 606-643 (2007).
2. Kaya, I.; Michno, W.; Brinet, D., et al. *Anal. Chem.* 89, 4685-4694 (2017).
3. Kaya, I.; Brinet, D.; Michno, W., et al. *ACS Chemical Neuroscience* 8, 2778-2790 (2017).
4. Aslund, A.; Sigurdson, C. J.; Klingstedt, T., et al. *ACS Chemical Biology* 4, 673-684 (2009).
5. Tycko, R. *Neuron* 86, 632-645 (2015).

TOr.05-3 VALIDATION STUDIES OF PHOSHOPEPTIDE ENRICHMENT PROCEDURES USING TiO₂-ATTACHED MEMBRANE FUNNEL-BASED SPRAY IONIZATION FOR IMAGING MASS SPECTROMETRY

Hung Y.-L. Winnie - Tak Wah Dominic Chan

The Chinese University of Hong Kong, Chemistry Department

Keywords: Imaging mass spectrometry, protein digestion, phosphopeptides, titanium dioxide, selective enrichment

Introduction

Imaging Mass Spectrometry (IMS) displays the distribution of different chemicals across a surface ^[1]. In this work, the analytical performance of an automated membrane funnel-based spray ion source was validated for phosphopeptides selective imaging mass spectrometry. Array of TiO₂-attached membrane funnels was fabricated for *in situ* clean-up and enrichment of phosphopeptides to alleviate matrix suppression effect in the ionization process(es) ^[2].

Methods

In-solution bovine beta-casein digest mixed with two phosphopeptide standards with different loadings were analyzed using funnel spray ionization mass spectrometry. In-funnel digestion of bovine beta-casein with different loadings was also performed. To evaluate the attainable spatial resolution of the IMS based on funnel spray ionization, "flower" figures were inkjet-printed onto Parafilm® membrane and images with different inter-funnel distances were obtained.

Results

The digestion derived phosphopeptides and the phosphopeptide standards were detected with good reproducibility. The analytical method exhibited reasonable linearity over a concentration range of 2.5 to 10 pmol for phosphoprotein loadings and 0.5 to 4.5 pmol for phosphopeptides loadings, demonstrating the feasibility TiO₂-attached membrane funnels in removing interfering non-phosphorylated peptides then selectively enriching phosphopeptides from complex protein digest samples. A mono-phosphorylated peptide was added as internal standard to the in-funnel digests before mass spectrometry analysis. This could reduce the experimental influence and improve the linearity of the analytical method. The best attainable spatial resolution for the IMS experiment was 200 microns and there was no detectable cross-contamination. Attempts to further reduce the inter-funnel distance using the current funnel fabrication method were hampered by the deformation of the membrane surface.

Conclusions

Phosphopeptides from protein digest were selectively detected with satisfactory detection limit, reproducibility and linearity using our developed workflow. IMS experiment with reasonable spatial resolution was demonstrated on this ion source. It is believed that this *in situ* single-pixel sample pretreatment method is a promising technique for the image studies of phosphoproteins of biological tissues.

Novel Aspect

An array of funnels could be used to pixelize a tissue surface for *in situ* single-pixel protein digestion and selective enrichment of phosphopeptides.

References

1. C. Wu, A. L. Dill, L. S. Eberlin, R. G. Cooks, D. R. Ifa, *Mass Spectrometry Reviews.*, 32, 218-243 (2013).
2. Li, W.; Chen, X. F.; Wong, Y. L. E.; Hung, Y. L. W.; Wang, Z.; Deng, L. L.; Chan, T. W. D. *Anal. Chim. Acta.*, 933, 97-102 (2016).

TOr.05-4 INFRARED LASER ABLATION MICROSAMPLING COUPLED WITH MALDI IMAGING

Kermit K. Murray⁽¹⁾ - **Fabrizio Donnarumma**⁽¹⁾ - **Kelin Wang**⁽¹⁾ - **Touradj Solouki**⁽²⁾ - **Michael Pettit**⁽²⁾

⁽¹⁾ Louisiana State University, Department of Chemistry, Baton Rouge - ⁽²⁾ Baylor University, Department of Chemistry, Waco

Keywords: MALDI, Imaging, Laser, Ablation

Introduction

Matrix-assisted laser desorption/ionization (MALDI) mass spectrometry imaging (MSI) is a powerful method for determining the location of biomolecules in tissue; however, protein identification and quantification remains challenging. The goal of this project is to develop an imaging workflow that combines MALDI imaging with laser ablation microsampling for liquid chromatography tandem mass spectrometry.

Methods

In the combined workflow, MALDI imaging is used to identify regions of interest (ROI) from intact proteins. The ROI are sampled using infrared laser ablation and the captured material is analyzed by LC-MS/MS using data independent acquisition to identify and quantify the proteins. The data are cross-correlated to identify the localized proteins in the MALDI images.

Results

Development of the combined approach is aimed at creating an automated system for ablation and capture and using it in a coupled workflow of MALDI imaging and LC MS/MS analysis. The infrared laser ablation and capture system uses a mid-infrared optical parametric oscillator laser with a custom reflective objective that has a large working distance and good numerical aperture. We have developed custom positioning software that allows MALDI MSI heat maps to be overlaid on camera images to co-register ROI ablation with the IR laser. Tissue sections are mounted on conductive microscope slides and either consecutive sections or MALDI analyzed sections can be used. Laser ablated proteins are digested with magnetic capture beads and the peptides released for analysis with a Waters nanoAcquity UPLC system coupled to a Synapt G2-HDMS.

Conclusions

MALDI MSI coupled with region-specific laser ablation sampling for LC MS/MS is a fast and versatile approach for spatially resolved tissue proteomics. We have demonstrated that proteins can be identified from spatially localized regions and are developing new methods for correlating the intact proteins observed in MALDI with the proteins identified by tandem mass spectrometry.

Novel Aspect

Coupled MALDI imaging with high precision infrared laser ablation capture for LC MS/MS for protein identification and quantification.

TO.r.05-5 PULSED COLD PLASMA FOR POST-IONIZATION IN MALDI-MS IMAGING

Jens Soltwisch⁽¹⁾ - Ulrich Röhling⁽²⁾ - Klaus Dreisewerd⁽¹⁾

⁽¹⁾ University of Münster, Institute for Hygiene, Münster - ⁽²⁾ University of Münster, Institute for Medical Physics and Biophysics, Münster

Keywords: Plasma-Postionization, MALDI-MSI, Enhanced sensitivity

Introduction

In MALDI, only a minor fraction of desorbed molecules is concomitantly ionized and, therefore, available for mass spectrometric analysis. To mitigate this inherent shortcoming, we here present a new strategy that utilizes a gas discharge for post ionization (PI) inside an intermediate pressure ion source. The RF-driven dielectric barrier discharge (DBD) is pulsed and timed to directly interact with the MALDI plume, leading to effective PI of matrix and analyte.

Methods

The DBD unit was integrated into the sample holder assembly of a SpectroGlyph dual funnel MALDI/ESI ion source, coupled to an Orbitrap MS. Consisting of an insulated antenna and a metal counter electrode, the DBD produces a spatially confined plasma right above the surface of the MALDI sample. Custom-made electronics were used to drive and trigger the discharge. MALDI-DBD-MS imaging experiments were performed with different matrices and in both polarities.

Results

In the first part of the presentation, concept and considerations for the design of the plasma region and RF-electronics for pulsed DBD-PI inside the intermediate pressure ion source will be described. This is followed by results from fundamental studies on the influence of relevant input parameters like DBD pulse duration, peak-to-peak voltage, delay between laser and DBD pulses, and the effect of gas pressure and its composition on the plasma. Careful optimization of laser and plasma parameters resulted in PI of phospho- and glycolipids with low degree of fragmentation. Under these optimized conditions, the DBD-MALDI mass spectra largely resemble those as produced by laser-induced PI in a MALDI-2 setting [1,2], i.e., a boost of protonated/deprotonated lipids and further metabolites is obtained. First DBD-MALDI imaging results obtained with mouse brain sections will be presented.

Conclusions

Pulsed DBD-PI holds great potential for increasing the sensitivity of MALDI-MSI in intermediate pressure ion sources at low costs. Synchronized to laser-generated material ejection, short bursts of DBD, suitably tailored in length and intensity, lead to a strong post ionization of matrix molecules and subsequent charge transfer reactions to neutral analytes within the evolving MALDI-plume.

Novel Aspect

A low-cost pulsed dielectric barrier discharge assembly for effective post ionization of lipids in an elevated pressure MALDI MS ion source

References

1. J. Soltwisch et al., *Science*, 348, 211-214 (2015)
2. S. R. Ellis et al. *Chem. Commun.*, 53, 7246-7249 (2017)

TOr.06 - Post Translational Modifications

Chairs: Joseph A. Loo, Fulvio Magni

TOr.06-1 **Keynote: RAPID, SENSITIVE AND QUANTITATIVE PHOSPHOPROTEOMICS FOR LARGE-SCALE CELL PERTURBATION ANALYSIS**

Jesper Velgaard Olsen

University of Copenhagen, Faculty of Health and Medical Sciences

Keywords: Phosphoproteomics, signaling, kinase, DIA, Ti-IMAC

Introduction

Protein activity is mainly modulated by dynamic reversible post-translational modifications (PTMs) such as site-specific phosphorylation, which regulates essentially all cellular processes. Despite continuous improvements, global analysis of protein phosphorylation is still challenging due to its sub-stoichiometric nature and low abundance.

Methods

We present a high-throughput IMAC based enrichment method to routinely quantify ~7000 phosphopeptides in 15 min of LC-MS/MS time on a Q Exactive HF-X MS [1]. But we have now reached the sequencing speed limit for data-dependent analysis (DDA) for phosphoproteomics. To overcome this, we have used data-independent analysis (DIA) for phosphoproteomics based on project-specific spectral libraries.

Results

To directly compare DDA and DIA in a biological setting, we performed a global kinase-substrate analysis, where we treated EGF stimulated RPE1 cells with three different MEK inhibitors. Phosphopeptides were enriched from only 200ug of peptides per sample using Ti-IMAC in 96-well format and analyzed with both DDA and DIA using label-free quantification. For the DIA analysis we generated a project specific spectral library containing ~55,000 unique phosphopeptides. From the DDA analysis we quantified ~7000 phosphopeptides per sample, whereas from the DIA analysis, we were able to quantify ~20,000 phosphopeptides per sample with few missing values. From the DIA runs we were also able to reproducibly identify ~400 phosphotyrosine-containing peptides in each sample, which is comparable to antibody-specific phosphotyrosine enrichments. Reassuringly, we found the same biological kinase signature using both DDA and DIA with strong overrepresentation of proline-directed ERK substrate motif sites among the inhibited phosphorylation sites.

Conclusions

DIA holds great promise for rapid phosphoproteomics profiling with deep coverage and accurate label-free quantification with few missing values.

Novel Aspect

Rapid DIA-based quantitative phosphoproteomics of hundreds of cell states

References

1. Kelstrup C.D., Bekker-Jensen D.B., Arrey T.N., Hoglebe A., Harder A., Olsen J.V. J. *Proteome Res.* 2018 17(1):727-738. doi: 10.1021/acs.jproteome.7b00602.

TOr.06-2 COMPREHENSIVE NISTmAb REFERENCE SPECTRAL LIBRARIES FOR RELIABLE AND ROBUST IDENTIFICATION OF POST-TRANSLATIONAL MODIFICATIONS IN HUMAN AND RECOMBINANT IGG ANTIBODIES

Qian Dong - Yuxue Liang - Xinjian Yan - Sanford Markey - Yuri Mirokhin - Dmitrii Tchekhovskoi - Tallat Bukhari - Stephen E. Stein

Biomolecular Measurement Division, National Institute of Standards and Technology, 100 Bureau Drive, Stop 8362, Gaithersburg, Maryland 20899, United States

Keywords: monoclonal antibodies, post-translational modifications, spectral library, 1D/2D LC-MS, glycation

Introduction

We report an alternative peptide and PTM identification method based on matching experimentally derived spectra in a comprehensive library of the NISTmAb IgG1 κ . We show that this method is advantageous over the database tools for reliably determining low levels of biologically modified species in antibody drugs. The simultaneous detection of these modifications, challenging for conventional methods, can be achieved without using special settings or procedures.

Methods

500 μ g of NISTmAb was denatured using various chemicals and solvents in 50 μ L of 100 mmol/L Tris HCl buffer at room or 85 °C. Parameters for 13 protocols were used during denaturing, reduction, and alkylation of the NISTmAb. Each sample was digested with trypsin for 0.25, 2, and 18h. 1D/2D LC-MS was used for various digestion conditions. Initial identifications were made using MS-GF+ search engine against a fasta-file containing the sequence.

Results

We created a high-resolution mass spectral library of all identified peptides produced in the tryptic digests of a humanized IgG1 κ mAb reference material (NISTmAb). The reference library contains over 12,600 high-quality tandem spectra of more than 3,300 peptide ions identified and validated by accurate mass, differential elution pattern, and expected peptide classes in peptide map experiments. These include a variety of biologically modified peptide spectra and method-induced artifacts obtained from 1D/2D LC-MS/MS analyses. The 12 varieties of spectra included in the library are derived from peptides of unmodified major classes, N-linked glycosylated, oxidized, deamidated, modified N-/C terminal, glycated, metal-ion adducted, overalkylated, and other types of modifications. A complete glycation profile was obtained for the NISTmAb with spectra for 58% and 100% of all possible glycation sites in the heavy and light chains, respectively. The site-specific quantification of methionine oxidation in the protein is described.

Conclusions

The library of annotated spectra of all NISTmAb peptide ions acquired in this work, along with the library search software, is freely available. The general utility of this library is demonstrated by the analysis of Humira, where 60% sequence coverage for constant regions of this commercial mAb. The NISTmAb library platform may be used as a tool for facile identification of the primary sequence and PTMs for human and recombinant IgG antibodies.

Novel Aspect

Facile, reliable, and simultaneous identification of low levels of biological modifications in mAbs by advanced NISTmAb reference material spectral library.

TOr.06-3 SEAGAL: THE METHOD TO FISH FOR UNCOMMON PHOSPHORYLATIONS

Vladimir Gorshkov⁽¹⁾ - Yaseelan Palarasah⁽²⁾ - Karsten Skjødt⁽²⁾ - Frank Kjeldsen⁽¹⁾

⁽¹⁾ Department of Biochemistry and Molecular Biology, University of Southern Denmark, Odense -
⁽²⁾ Department of Cancer and Inflammation Research, University of Southern Denmark, Odense

Keywords: phosphorylation, enrichment, gallium, immunopurification

Introduction

Traditional methods for phosphopeptide purification are reported to introduce strong bias and hence detected phosphorylations only partially reflect the true phosphoproteome distribution of the cell [1]. Importantly, part of phosphoproteome, especially, N-linked phosphorylations (*i.e.* histidine, arginine, lysine) are generally out of the scope of conventional enrichment methods, despite high interest in those type of phosphorylation.

Methods

The mixture of phosphopeptides was modified by the recently developed by dimetallic gallium complex [2] in 20 mM triethylammonium bicarbonate (pH 8). Modified peptides were enriched by immunopurification using an immobilized antibody against the introduced metal complex and analyzed by LC-MS/MS. Raw data analysis was performed in Thermo Proteome Discoverer 2.2, including home-built data deconvolution and pre-processing step and conventional database search.

Results

The developed dimetallic complex is highly specific for the phosphate moiety, irrespective of its nature, thus, the same method allows simultaneous enrichment of peptides having O-linked (Ser, Thr, Tyr) and N-linked (His, Lys, Arg) phosphorylation. The modification and enrichment are performed at close to neutral pH (optimal range 5 – 8), thus, avoiding destabilization of N-linked phosphorylation and providing less bias of the phosphoproteome. Modified peptide and all fragments containing the modification display distinct isotopic pattern preventing efficient data analysis by conventional proteomics tools. To solve this issue, data pre-processing software was developed as a processing node for Thermo Proteome Discoverer. The source code and complete description of algorithms are made available for the community. Alternatively, the introduced modification could be removed by pyrophosphate in a simple one-step process compatible with conventional LC-MS-based proteomics.

Conclusions

The proposed phosphopeptide enrichment approach is based on specific phosphate-capturing dimetallic complex and immunopurification of complexed phosphopeptides. The enrichment can be performed at neutral pH, preserving labile modifications. The present proof-of-concept study clearly indicates the potential of the method, however, further optimization is necessary to apply the method in biologically-relevant cases.

Novel Aspect

The presented method allows for obtaining comprehensive information on O- and N-linked phosphoproteome simultaneously at neutral pH.

References

1. Pan, C., Gnad, F., Olsen, J. V., Mann, M. *Proteomics*, 8, 4534-4546 (2008)
2. Svane, S., Kryuchkov, F., Lennartsson, A., McKenzie, C. J., Kjeldsen, F. *Ang. Chem. Int. Ed.*, 51, 3216-3219 (2012)

TOr.06-4 PROTEOME-WIDE EFFECTS OF SINGLET OXYGEN PRODUCED BY NEXT GENERATION IRIIDIUM ANTI-CANCER METALLODRUGS

Cookson K. C. Chiu - Pingyu Zhang - Pui Yiu Yuko Lam - Christopher A. Wootton - Mark P. Barrow - Peter J. Sadler - Peter B. O'Connor

Department of Chemistry, University of Warwick, Coventry

Keywords: anticancer, photodynamic therapy, proteomics, Fourier Transform Ion Cyclotron Resonance (FT-ICR) MS, post-translational modification

Introduction

New metal-based anticancer drugs have been the focus of next-generation chemotherapy. Novel Ir(III) based drugs have been designed to minimise side effects on healthy cells, while increasing selectivity and efficacy towards cancer cells by inducing reactive oxygen species (ROS) when triggered with blue visible light. The catalytic production of $^1\text{O}_2$ with Ir-complexes, and effects on target species within cells were investigated to uncover their mechanism of action on a proteome-wide scale. [1]

Methods

FT-ICR MS analysis was conducted using a 12T Bruker Solarix (Bruker Daltonik GmbH, Bremen, Germany) FT-ICR mass spectrometer. nLC-MS/MS experiments was conducted via coupling an EASY nano-LC II system (Proxeon, Hemel Hempstead, UK) to the Solarix FT-ICR MS, separation was achieved using a water/acetonitrile gradient and a C18 reverse phase column. All data was analysed with Data Analysis v4.2 (Bruker Daltonik GmbH, Bremen, Germany) and the Mascot database search engine.

Results

Activation of the Ir metallodrugs caused production of reactive singlet oxygen species and extensive oxidation of nearby biomolecules. A whole range of oxidation products, up to 6 per peptide, were observed. The high resolution and mass accuracy and extensive MS/MS capabilities of the FT-ICR MS allowed unambiguous identification of the oxidation sites and the modifications induced. $^1\text{O}_2$ was found to cause many methionine oxidation events (to sulfoxide and sulfone), but also oxidation of tryptophan residues, producing diagnostic kynurenine and 3-hydroxy kynurenine moieties. Histidine residues were oxidised to 2-oxo-histidine, which is commonly observed in $^1\text{O}_2$ oxidation events. Bottom-up nLC-MS/MS experiments of drug-treated cancer cells revealed two key targets; aldose reductase (AR) and heat shock protein 70 (Hsp70). Quantitative nLC-MS/MS showed the increase in oxidised products in the drug-treated samples, calculated as 3.0-fold up-regulation with AR and 5.8-fold up-regulation with Hsp70, together with 9 proteins up-regulated along the glycolysis pathway.

Conclusions

New Ir-based anticancer drug was designed to be highly and selectively cytotoxic to cancer cells when activated with visible blue light. Generation of $^1\text{O}_2$ was confirmed by observing diagnostic oxidised species in peptides. Damage to specific cellular proteins, AR and Hsp70 was discovered in cells, induced by photodynamic therapy. 9 proteins along the glycolysis pathway were found up-regulated, which happened in the mitochondria of cells, disrupting the energy source of cancer cells.

Novel Aspect

Top-down and bottom up FT-ICR MS/MS approaches shed light on new anticancer, catalytic photoactivated Iridium complex targets and proteome effects

References

1. Zhang P., Chiu C.K.C., Lam Y.P.Y., Wootton C.A., Barrow M.P., Sadler P.J., O'Connor P.B., *Angew. Chem. Int. Ed.*, 56, 14898-14902 (2017).

TOr.06-5 POLYAMINE MODIFICATIONS IN DIATOM BIOSILICA PROTEOMES

Alexander Milentyev⁽¹⁾ - Christoph Heintze⁽²⁾ - Nicole Poulsen⁽²⁾ - Nils Kroeger⁽²⁾ - Andrej Shevchenko⁽¹⁾

⁽¹⁾ Max Planck Institute of Molecular Cell Biology and Genetics, Shevchenko Lab, Dresden - ⁽²⁾ B-CUBE Center for Molecular Bioengineering, Kroeger Lab, Dresden

Keywords: diatoms, biosilica, post-translational modification, polyamines, AQC

Introduction

Diatoms are unicellular algae that use highly specialized proteins to produce nano silicified cell walls. These proteins, termed *silaffins*, share no homology across diatom species. Silaffins are heavily modified [1] and their lysine residues bear polyamine chains [2]. However, their structure and physiological role remain elusive.

Methods

Total biosilica extracts from *Thalassiosira pseudonana*, *T. oceanica* and *Cyclotella cryptica* species were isolated as described [3], subjected to acidic hydrolysis and derivatized with AccQ•Ultra Reagent (Waters). Lysine derivatives were quantified by LC-MS and confirmed by MS/MS on Q-Exactive HF mass spectrometer (Thermo Scientific). Modifications were mapped to protein sequences by GeLC-MS/MS [4] on OrbitrapVelos mass spectrometer (Thermo Scientific).

Results

We developed a method to identify and quantify polyamine-modified lysines and map them back to silaffin sequences. In total, 17 novel modifications were discovered, including acid-resistant phosphoester-containing polyamines. We demonstrated that the pattern of polyamine modifications reflects the phylogenetic proximity of the diatom species. Modified lysine residues were identified by polyamine-specific fragments in MS/MS spectra followed by iterative searches and deconvolution of raw MS/MS spectra. We localized 130 polyamine-modified sites in 26 proteins from three diatom species and revealed 3 consensus motives common to all three diatoms.

Conclusions

We report a novel approach to identify, quantify and localize 130 polyamine modifications in three diatoms and demonstrated that they occurred at conserved consensus motifs, despite full sequences of the modified proteins were not conserved.

Novel Aspect

Identification of consensus motifs of polyamine modifications

References

1. N. Kröger and N. Poulsen, *Annu. Rev. Genet.*, vol. 42, pp. 83–107 (2008)
2. M. Sumper, R. Hett, G. Lehmann, and S. Wenzl, *Angew. Chem. Int. Ed. Engl.*, 46, 8405–8 (2007)
3. N. Kröger, R. Deutzmann, C. Bergsdorf, and M. Sumper, *Proc. Natl. Acad. Sci. U. S. A.*, 97, 14133–8 (2000)
4. A. Shevchenko, H. Tomas, J. Havlis, J. V Olsen, and M. Mann, *Nat. Protoc.*, 1, 2856–60 (2006)

TOr.07 - Mass spectrometric analysis of forensic science evidence

Chairs: Simona Francese, Gloriano Moneti

TOr.07-1 **Keynote: EXTENDED MASS SPECTROMETRY WORKFLOW FOR THE ANALYSIS OF DRUGS OF ABUSE IN BIOLOGICAL FLUIDS AND TISSUES**

Gérard Hopfgartner

Life Sciences Mass Spectrometry, University of Geneva, Switzerland

Keywords: Drugs of Abuse, LC-MS, Data Independent Analysis, SWATH, Differential Ion Mobility Spectrometry, ExD

Introduction

Screening of Drugs of Abuse (DoA) in biological matrices is mainly based on triple quadrupole or high resolution mass spectrometers using data dependent acquisition workflows (DDA). The development of data independent acquisition such as SWATH, differential ion mobility (DMS) and electron induced dissociation (ExD) opens new possibilities for comprehensive and multidimensional for simultaneous qualitative and quantitative analysis.

Methods

Data were acquired on either a triple quadrupole linear ion trap (QTRAP5500, Sciex) equipped with DMS (Selexion, Sciex) or a quadrupole time of flight equipped with DMS and a prototype chimeric collision cell allowing to perform collision induced dissociation (CID) or electron induced dissociation (EID). Human urine samples were analyzed using reversed phase liquid chromatography.

Results

Various workflows have been developed for the direct analysis of DoA on tissues or in biological fluids such as urine using multidimensional MS approaches. For direct MS analysis using ambient ionization or for LC-MS modifier assisted DMS offers an additional separation dimension. Contrarily to pharmaceutical investigation toxicological screening cannot rely on the origin/nature of the drug neither on the intake time nor the frequency. Comprehensive LC-SWATH/MS allows to collect in the same analysis precursor ions as well as fragment ions. Information can be obtained either hypothesis driven or with the use of multivariate analysis post-acquisition. This provides a powerful tool for the screening of DoA as well as the monitoring of the exposome or to monitor the change in endogenous metabolites in human urine samples. In addition to CID electron induced dissociation of singly charged precursor ions (protonated or adducts) generates additional information for structural identification of metabolites.

Conclusions

The complexity of biological samples and the problematic of DoA call for improved methods for screening and quantification. Novel development in mass spectrometry such as DIA, ion mobility and EID with or without chromatography separation can be used to extend the toxicological knowledge as the intake of a DoA can also be monitored based on their contaminants.

Novel Aspect

Integration of data independent acquisition, ion mobility and electron induced dissociation in drug abuse screening.

References

1. Bonner R., Hopfgartner G., *Bioanalysis*. 8, 1735-50 (2016)

2. Baba T., Campbell J.L., Le Blanc Y.J.C, Hager J.W., Thomson B.A., *Anal. Chem.* 87, 1, 785-792 (2015)

TO.r.07-2 THE DEATH OF THOMAS CHATTERTON - AN LC-MSⁿ INVESTIGATION

Paul Gates

University of Bristol, School of Chemistry, Bristol

Keywords: morphine, degradation, forensics, LC-MS/MS

Introduction

Thomas Chatterton (1752 – 1770) was a literary genius who tragically died in unexplained circumstances [1]. At the time, the Coroner claimed that he “swallowed arsenic in water and died, in consequence thereof” despite no arsenic being found at the scene. In his belongs was a memorandum book, which had a dark stain running through its pages. A police forensics report in 1947 stated that the stain “gives a positive reaction for opium alkaloids” but gives no further details.

Methods

Herein, is a study to try to shed some light on Chatterton’s death using high resolution LC-MSⁿ on an Orbitrap Elite instrument. Preliminary tests were performed with codeine (an opium alkaloid) to test if it was possible to recover it from paper after being deposited. This test proved successful, so then a sample of a stained page from the memorandum book was extracted using methanol and dichloromethane and analysed by the same methods.

Results

The preliminary study showed the recovery of codeine from office paper. The extracts also contained unidentifiable components probably from the paper itself, as well as the common phthalate and siloxane contaminants. Codeine was recovered intact as well some possible degradation products, but these were not investigated further. In the main study, a 2 mm² sample of a stained page from the memorandum book was extracted and analysed by both direct infusion ESI-MS and HPLC-MSⁿ. As expected, from a sample of this age, the extract was very complex, and no signal was observed for any intact opium alkaloids. The LC-MS analysis, however, showed a cluster of at least 15 peaks with molecular weights in the 200-300 *m/z* range. 8 of these have been confirmed by accurate mass MS/MS analysis to be known degradation products of morphine [2]. A further 7 components have also been identified as being possible unpublished degradation products of morphine, codeine or thebaine.

Conclusions

This study demonstrates the extraction and identification of at least 8 known morphine degradation products from a 250-year-old book by high-resolution HPLC-MSⁿ. The presence of a further 7 opium alkaloid degradation products is also proposed. These results indicate that the stain originally contained opium alkaloids, most probably in the form of laudanum – a popularly prescribed medical treatment of the day.

Novel Aspect

The demonstration of modern analytical methods to the study of a ‘mystery’ stain in a 250-year-old book found at the scene of an unexplained death in 1770.

References

1. Gregory G., Kippis A. (ed), *Biographia Britannica*, 4, 573-619 (1789)
2. Rapoport H., *Journal of Organic Chemistry*, 13, 714-721 (1948)

TOr.07-3 CLASS IDENTIFICATION OF NOVEL PSYCHOACTIVE SUBSTANCES USING MASS SPECTRAL LIBRARY SEARCHING WITH HIERARCHICAL CLUSTERING

Arun Moorthy⁽¹⁾ - **William Wallace**⁽¹⁾ - **Anthony Kearsley**⁽²⁾ - **Dmitrii Tchekhovskoi**⁽¹⁾ – **Stephen E. Stein**⁽¹⁾

⁽¹⁾ National Institute of Standards and Technology, Mass Spectrometry Data Center, Gaithersburg - ⁽²⁾ National Institute of Standards and Technology, Applied and Computational Mathematics Division, Gaithersburg

Keywords: Designer Drug, Forensics, Hierarchical Clustering, Library Searching, Reference Data

Introduction

Novel psychoactive substances have introduced significant identification challenges to the forensic science community. Developing methods that can correctly categorize emerging drug analogs is a pressing need. This presentation describes a procedure for obtaining a class-identification *map* of an unknown compound. The resulting classification map can be used to compute the uncertainty of the class-identification.

Methods

The *hybrid search* is a new library search method that considers fragment-ions along with neutral-losses when computing spectral similarity [1]. A query spectrum is first searched against a library to generate a preliminary hit-list. The *hybrid* dot-product calculation is used to generate a dissimilarity matrix from the spectra in the preliminary hit-list, which is then clustered using our novel multi-objective agglomerative-hierarchical clustering algorithm.

Results

The result of this mass spectral library search and clustering procedure is a dendrogram from which the likely classification of the query compound can be proposed. This dendrogram can be used to compute the uncertainty of the proposed class identification. As a numerical test, twenty-two SWGDRUG Library [2] compounds were assembled into a sub-library. It contained spectra for five analogs each of JWH-167, fentanyl, MDMA, and cathinone, and two additional spectra of “unknown” compounds. A preliminary implementation of the multi-objective classification method was able to classify the first unknown as a fentanyl analog with a certainty of 56%, and as a fentanyl, JWH-167, or cathinone analog with 93% certainty. The second unknown was only classifiable as a fentanyl, MDMA, or cathinone analog, but with 99% certainty.

Conclusions

This presentation describes a procedure for obtaining a class-identification for an unknown seized compound, using the new *hybrid search* to compute inter-spectral similarity, and a novel multi-objective scheme for agglomerative hierarchical clustering. The resulting classification tree can be used to compute the uncertainty of the class-identification.

Novel Aspect

The hybrid dot-product for spectral similarity, a new multi-objective hierarchical clustering algorithm, and computation of uncertainty for class identification.

References

1. Moorthy A.S., Wallace W.E., Kearsley A.J., Tchekhovskoi D.V., Stein S.E., *Analytical Chemistry*, 89, 13261-13268 (2017)
2. www.swgdrug.org/ms.htm

For information please contact: william.wallace@nist.gov

TOr.07-4 IDENTIFYING WOOD FROM PROTECTED SPECIES WITH DART FINGERPRINT SPECTRA AND THE SOFTWARE PLATFORM OMFF

Sandra Martínez-Jarquín⁽¹⁾ - **Cady Lancaster**⁽²⁾ - **Ed Espinoza**⁽³⁾ - **Robert Winkler**⁽¹⁾

⁽¹⁾ CINVESTAV, Biochemistry and Biotechnology, Irapuato - ⁽²⁾ U.S. Forest Service, National Fish and Wildlife Forensic Laboratory, Utah - ⁽³⁾ U.S. Forest Service, National Fish and Wildlife Forensic Laboratory, Oregon

Keywords: DART-MS, wood identification, software, fingerprinting, ambient ionization

Introduction

According to the Food and Agricultural Organization (FAO), global production of major wood products increases at a 3-6% per year. Some commercially valuable wood treated products are illegally used from endangered tree species, that are not easily identified. DART-MS has been used to differentiate morphologically similar species with great results[1-3]. We want to build a software with a scalable reference database for forensic applications.

Methods

Reference spectra are generated from wood libraries at the Forensic Laboratory for Wildlife and Fish (US Forest Service) from different tree species and samples. An SQL reference database was constructed with 4,950 spectra from 301 different species. The Open Mass Fingerprinting Framework (OMFF) was designed as a webpage for easy distribution, with html and bootstrap at the front end and Django to improve security of the back end. Multivariate statistical analysis was implemented for data processing and identifying specific chemical markers.

Results

A software for tree identification based on the DART-MS chemical profile data was constructed (OMFF). DART-MS provides a fast analysis which requires minimal time and amount of sample making it ideal for processed wood products (guitars, furniture, etc.). The program is user-friendly, open source webpage that includes a reference database for identifying of wood samples. Specific chemical markers are obtained from DART-MS fingerprint reference spectra and compared with unknown samples. The results provide a list of species with a comparison probability score. The software produces statistically meaningful information. While the number of samples for some species in the reference database is small, the system can be expanded. The program was designed for non-expert users; therefore, it can be implemented at critical points such as customs and ports.

Conclusions

The implementation of the OMFF software will help diminish traffic of protected wood species. CITES provides tools for wood identification such as morphological, DNA and microscopy characterization. Using all the available tools will improve identification confidence levels for wood species.

Novel Aspect

OMFF provides a reference database and reliable classification mechanisms for non-expert users to identify morphologically similar wood samples.

References

1. McClure, P.J., Chavarria, G.D., & Espinoza, E., *Rapid Commun. Mass Spectrom.*, 29, 783–788 (2015)
2. Evans, P.D., Mundo, I.A., Wiemann, M.C., Chavarria, G.D., McClure, P.J., Voin, D., & Espinoza, E.O., *IAWA J.*, 38, 266–281 (2017)
3. Deklerck, V., Finch, K., Gasson, P., Van den Bulcke, J., Van Acker, J., Beeckman, H., & Espinoza, E., *Rapid Commun. Mass Spectrom.*, 31, 1582–1588 (2017)

TO.r.07-5 GREEN FORENSICS - ISOTOPIC TRACEABILITY IN HORTICULTURE. GEOGRAPHICAL ORIGIN OF ROSES

Eva de Rijke⁽¹⁾ - Meyci Trisna⁽²⁾ - Suzan Verdegaal⁽³⁾ - Runa Magnússon⁽¹⁾ - Jorien Schoorl⁽¹⁾ - Chiara Cerli⁽¹⁾ - Chris de Koster⁽²⁾

⁽¹⁾ University of Amsterdam, Institute for Biodiversity and Ecosystem Dynamics, Amsterdam - ⁽²⁾ University of Amsterdam, Swammerdam Institute of Life Sciences, Amsterdam - ⁽³⁾ VU University Amsterdam, Faculty of Earth and Life Sciences, Amsterdam

Keywords: forensics, isotope ratio mass spectrometry, GC-MS, roses, geographical origin

Introduction

Roses are the most important Dutch cut flowers and are grown by companies under a global GAP (Good Agricultural Practice) certificate that is used by e.g. European supermarket organizations. Certified companies have strict requirements for safe, sustainably grown, high-quality and traceable products. Next to the Netherlands, Kenya and Ecuador are the main exporting countries for roses. The question rises whether roses of unknown origin, which do not meet these strict requirements, or roses from areas with a phytosanitary export ban due to occurrence of certain plant diseases, are exported illegally under this certificate to the EU. Therefore, UvA, Naktuinbouw and De Ruiter BV develop isotope analysis technology for provenance determination of roses. The stable isotopic composition of H, C, N and O from plant materials depends on regional geological and climatic patterns, and therefore provides a regional fingerprint of the plant.

Methods

In this study, we focus on development of methods for isotope ratio and gas chromatography (GC) mass spectrometry analysis of stems and water from roses that were harvested monthly throughout the 9-month season. Elemental analysis (EA) coupled to isotope ratio mass spectrometry (IRMS) was applied to investigate the use of $\delta^{13}\text{C}$, $\delta^{15}\text{N}$ and total C and N profiles in the stems as a fingerprint for the country of origin. Furthermore, the $\delta^{18}\text{O}$ isotopic composition of the xylem water in the rose stems, extracted via freeze-drying, was linked to $\delta^{18}\text{O}$ of irrigation water and the GNIP (Global Network of Isotopes in Precipitation) IAEA database, to determine differences between roses from different geographical origins. A third method used was quantitative GC-MS analysis and compound specific GC-IRMS analysis of leaf wax alkanes, extracted via solid phase extraction.

Results

The best separation of the groups was obtained from the water analysis using a combination of the $\delta^2\text{H}$ and $\delta^{18}\text{O}$ values of the rose and irrigation source water. The optimal time of freeze drying was at least 22 hours. $\delta^2\text{H}$ and $\delta^{18}\text{O}$ values of source water do not vary much over the season and the variation within one batch is low, furthermore they correspond well to the global mean water line (GMWL). The results show a correlation between the source water and the local long-term GNIP averages (Pearson's R = 0.9168 for $\delta^2\text{H}$ and 0.9029 for $\delta^{18}\text{O}$) and between the stem- and source water (Pearson's R = 0.9959 for $\delta^2\text{H}$ and 0.9921 for $\delta^{18}\text{O}$), with constant isotope fractionation between the source and the rose water data of about $-4.\delta^{13}\text{C}$ and $\delta^{15}\text{N}$ values of bulk rose stem analysis with EA-IRMS can also be used to distinguish roses from different origins, especially Dutch roses which were clearly separated from Ecuadorian and Kenyan samples. The measured $\delta^{15}\text{N}$ isotope values also correspond well to the reported literature values [1], even though other plants were used in this study. GC-MS analysis of the relative concentrations of C_{27} , C_{29} , C_{31} and C_{33} in leaf wax extracts showed clustering according to the country of origin. In a previous study, similar results were obtained for the same alkanes in bell pepper [2, 3].

Conclusions

Three analytical strategies were investigated, based on IRMS analysis of $\delta^2\text{H}$ and $\delta^{18}\text{O}$ of rose and

irrigation source water, $\delta^{13}\text{C}$ and $\delta^{15}\text{N}$ of rose stems, and on quantitative GC-MS analysis of leaf wax *n*-alkanes. All strategies were able to discriminate between Dutch, Ecuadorian and Kenyan roses that were monitored monthly throughout the season. The source water and long-term regional GNIP data showed good correlation with the rose water data, with constant isotope fractionation of about -4. These multi-isotopic methods are therefore suited as a tool to verify the geographical origin of roses in order to prevent fraud and phytosanitary problems.

Novel Aspect

This is the first time a combined strategy of multi-isotopic EA-IRMS and quantitative GC-MS methods has been developed as a tool to verify the geographical origin of roses.

References

1. Amundson R., Austin A. T., Schuur E. A. G., Yoo K., Matzek V., Kendall C., A. Uebersax, D. Brenner, and W. T. Baisden, *Global Biogeochemical Cycles*, 17, 1031 (2003).
2. de Rijke E., Schoorl J., C. Cerli, Vonhof H., Verdegaal S., Vivo-Truyols G., Lopatka M., Dekter R., Bakker D., Ebskamp M., de Koster C.G., *Food Chemistry*, 204, 122 (2016).
3. de Rijke E., Fellner C., Westerveld J., Lopatka M., Cerli C., Kalbitz K., de Koster C.G., *Analytical and Bioanalytical Chemistry*, 407 5729(2015).

TOr.08 - Ion activation and dissociations

Chairs: Vicki Wysocki, Gianluca Giorgi

TOr.08-1 UNDERSTANDING FRAGMENTATION OF COMPLEX IONS USING SURFACE-INDUCED DISSOCIATION EXPERIMENTS

Julia Laskin

Purdue University, Dept. of Chemistry, West Lafayette, IN 47907

Keywords: surface-induced dissociation, threshold energy, activation entropy, peptide, metal cluster

Introduction

Surface-induced dissociation (SID) is a powerful tool both for analytical applications and for studies of gas-phase fragmentation of complex ions. We have developed experimental approaches for studying the energetics, dynamics, and kinetics of large ion dissociation using SID, which in combination with comprehensive modeling enable mechanistic understanding of complex ion fragmentation difficult to obtain using other experimental techniques.

Methods

Experiments were performed on a specially designed Fourier transform ion cyclotron resonance mass spectrometer (FT-ICR) configured for time- and collision energy-resolved SID experiments. The experimental results were modeled using RRKM and the analytical form of the internal energy deposition function that adequately describes energy transfer in collisions.

Results

SID experiments were combined with resonant ejection of selected fragment ions. The approach was first examined using well-characterized peptide systems and subsequently used to explore the kinetics of the gas-phase fragmentation of ligated metal clusters. When combined with RRKM modeling, these experiments facilitate detailed understanding of the competition between different fragmentation pathways.

Conclusions

SID experiments highlight an important role of entropy in dissociation of complex ions. Resonant ejection experiments identify different populations of stable structures and inform on the kinetics of the competing fragmentation reactions.

Novel Aspect

SID provided unique insights into the energetics of complex ion fragmentation.

TO.r.08-2 INVESTIGATION OF HEMICRYPTOPHANE HOST-GUEST BINDING ENERGIES USING ADVANCED TANDEM MASS SPECTROMETRY

Parisa Bayat⁽¹⁾ - David Gatineau⁽²⁾ - Denis Lesage⁽¹⁾ - Vincent Robert⁽³⁾ - Alexandre Martinez⁽⁴⁾ - Richard B. Cole⁽¹⁾

⁽¹⁾ Sorbonne Université, IPCM, Paris - ⁽²⁾ Univ Grenoble Alpes, DCM, Grenoble - ⁽³⁾ Université de Strasbourg, LCQ, Strasbourg - ⁽⁴⁾ Aix Marseille Univ, Equipe Chirosciences, Marseille

Keywords: Host-Guest Chemistry; cryptophanes; blackbody infrared radiative dissociation; collision induced dissociation; MassKinetics

Introduction

In advancing host-guest (H-G) chemistry, considerable effort is spent to synthesize host molecules with specific properties including selectivity and adjustable affinity. An important step in the process is the characterization of binding strengths of the H-G complexes. Here, we present a mass spectrometry based multimodal approach to obtain bond dissociation energies (BDEs) of two hemicryptophane cages with three biologically-relevant guest molecules.

Methods

Two model hemicryptophane cages have been selected as the host molecules and their binding properties towards three biologically important guests: choline (C), acetylcholine (AC) and betaine (B) have been studied in the gas phase. Quantitative comparison of the binding dissociation energies of different H-G pairs was achieved using a combination of BIRD and high-pressure CID mass spectrometric techniques with Rice–Ramsperger–Kassel–Marcus (RRKM) modeling.

Results

Obtained BDE values fall within a rather narrow range for all host-guest complexes. Among the informative trends and differences extracted from the data, we note that complexes of cage 1 exhibit higher BDE's, and thus are all more stable than those of cage 2. This trend establishes that the linkers (naphthyl for cage 1 and phenyl for cage 2) play important roles in determining the strength of the H-G complexes. Another interesting comparison is the relative stability ranking between various guests for a given host cage. Notably, for both cages 1 and 2, the order of stability of H-G complexes is as follows: choline > acetylcholine > betaine. DFT calculations were performed to gain information concerning possible binding modes and relative stabilities of different non-covalent H-G interactions. The trend and the CTV of the hemicryptophane adopted an approximately C₃ symmetry. The arrangement of the phenyl or naphthyl linkers for the six H-G complexes gives compact structures that allow the insertion of the ammonium moiety of the guest inside the cage.

Conclusions

The multimodal approach described here, using the BIRD technique in combination with high pressure-CID breakdown curves and RRKM modeling, is useful not only for estimation of bond dissociation energies, but also for accurate surveying of the influence of different functional groups on H-G binding properties. Utilization of this approach can provide a wealth of

information pertinent to chemists who seek to design new H-G systems with well-honed properties.

Novel Aspect

Two hemicyptophane host molecules are assessed for their binding properties towards three biologically important guests in the gas phase.

TOr.08-3 TANDEM TOF MASS SPECTROMETRY OF BIOMOLECULES: FROM DUAL STAGE TO CURVED FIELD AND PARABOLIC ANALYZERS AS SECOND STAGE

Ernst Pittenauer - Günter Allmaier

TU Wien, TU Wien/Institute of Chemical Technologies and Analytics, Vienna

Keywords: MALDI, tandem TOF instrumentation, high-energy CID, charge-remote fragmentation

Introduction

By the end of the millennium high-energy CID were moved from FAB/LSIMS tandem sector to MALDI tandem TOF devices. Surprisingly, this generation of instruments exhibits extremely different instrumental features (lasers, collision energies, reflectrons and flight path lengths). In order to demonstrate similarities/differences among high-energy CID spectral data from peptides, lipids and oligosaccharides it is of great interest for future applications [1-4].

Methods

MALDI high-energy CID spectra of selected peptides, triacylglycerols, glycerophospholipids and oligosaccharides were acquired on an AB Sciex 4700/4800 ($E_{\text{LAB}}=1$ keV), a Bruker Daltonics ultrafleXtreme ($E_{\text{LAB}}=8$ keV), a Shimadzu Axima TOF² ($E_{\text{LAB}}=20$ keV), a Shimadzu MALDI 7090 ($E_{\text{LAB}}=20$ keV) and a JEOL SpiralTOF/TOF instrument ($E_{\text{LAB}}=20$ keV). MALDI matrices used included α -cyano-4-hydroxy cinnamic acid (only peptides) or 2,4,6-trihydroxy acetophenone.

Results

Whereas tandem TOF instruments fitted with a limited energy acceptance dual stage reflectron (AB Sciex, Bruker) are typically operated with a rather low collision energy (1-8 keV), the ones with a wide energy acceptance curved field or offset parabolic reflectron (Shimadzu, Jeol) allow a collision energy of 20 keV. Protonated peptide molecules only exhibit side chain charge-remote fragmentations (d-, v-, w-ions) at 20 keV collisions allowing distinction between isomeric leucine and isoleucine. Triacylglycerols and glycerophospholipids show - depending on the type of precursor ion (protonated, sodiated) - diagnostic charge-remote fragmentation of the alkyl chains best seen at 20 keV collisions. For sodiated triacylglycerols a clear differentiation between *sn1-/sn3-* versus *sn2-* substituents could be made, hardly seen at 1-8 keV collisions. In contrast, sodiated oligosaccharides do not seem to require such high collision energies as there is virtually no difference in the fragmentation pattern in the CID spectra taken at 1, 8 or 20 keV.

Conclusions

Only instruments fitted with wide energy acceptance reflectrons (curved field and offset paraboloid reflectron) yield true high-energy CID spectra of biomolecules with typical product ions in high abundance (peptide side-chain fragmentation, lipid charge-remote fragmentation). Still a major issue for the observation of unwanted PSD ions in most tandem TOF instruments is a long flight path (~4 m) obscuring CID ions due to their lower intensity (a factor of 10-20).

Novel Aspect

For the first time different designs of MALDI tandem TOF devices fitted with high-energy CID cells are compared for their usefulness in structural elucidation of biomolecules.

References

1. Pittenauer E., Allmaier G., Comb Chem High Throughput Screen, 12, 137 (2009).
2. Pittenauer E., Allmaier G., J. Am. Soc. Mass Spectrom., 20, 1037 (2009).
3. Pittenauer E., Allmaier G., Int. J. Mass Spectrom., 301, 90 (2011).
4. Pittenauer E., Rehulka P., Winkler W., Allmaier G., Anal. Bioanal. Chem., 407, 5079 (2015).

TOr.08-4 **SPECIFIC RADICAL-DRIVEN DISSOCIATION OF HYDROGEN-EXCESS AND -DEFICIENT RADICAL PEPTIDES PRODUCED BY MALDI IN-SOURCE DECAY MASS SPECTROMETRY**

Mitsuo Takayama

Yokohama City University, Graduate School of Nanobioscience, Yokohama

Keywords: Radical-Driven Dissociation, MALDI-ISD, Hydrogen-Excess/Deficient Peptides, Mechanism

Introduction

Radical-driven dissociation (RDD) using ultraviolet, electron transfer, high-energy collision and redox matrix provides us practical information for structure elucidation of organic chemicals, peptides and proteins. Using MALDI MS with redox matrix here we report detailed mechanisms of the formation of hydrogen-excess/deficient radical peptides, $[M+H]^{\cdot}$ and $[M-H]^{\cdot}$, and specific RDD reactions.

Methods

MALDI-ISD experiments were performed with a TOF MS AXIMA-CFR (Shimadzu, Kyoto, Japan). Redox matrices used were 2,5-dihydroxybenzoic acid (2,5-DHB) and 5-amino-1-naphthol (5,1-ANL) for H-donation and 5-nitrosalicylic acid (5-NSA) and 4-nitro-1-naphthol (4,1-NNL) for H-abstraction [2]. The analytes used were synthetic peptides RGV and EGV having Arg at N-terminus (RGVKVDYGPAMAFAPADLAAEF, RGV) and C-terminal side (EGVKVDYGPAMAFAPADLAAEF, EGV).

Results

ISD spectrum of RGV with H-donating matrix 5,1-ANL merely gave $[c+2H]^+$ ions originated from cleavage at N-Ca bond [1], while the use of 2,5-DHB resulted in $[c+2H]^+$ and $[a]^+$ ions and some $[b]^+$ and $[d]^+$ ions. The ISD spectrum of EGV with 5,1-ANL showed preferential $[z+2H]^+$ ions and minor $[y+2H]^+$ and $[c+2H]^+$ ions, while the spectrum with 2,5-DHB gave preference $[y+2H]^+$ ions and minor $[z+2H]^+$, $[x]^+$ and $[c+2H]^+$ ions. The use of 5,1-ANL resulted in only H-excess radical peptide ions $[M+2H]^{\cdot}$ having the radical (unpaired electron) located on carbonyl carbon of the backbone, and the radical ions $[M+2H]^{\cdot}$ can lead to $[c+2H]^+$ ion and radical intermediate ion $[z+H]^{\cdot}$, although the intermediate never be observed in ISD spectra. On the other hand, the use of H-abstracting matrix 4,1-NNL resulted in preferential $[a]^+$ ions and some $[b]^+$ and $[d]^+$ ions for RGV and preferential $[y]^+$ ion for EGV.

Conclusions

The use of naphthalene skeletal redox matrices 5,1-ANL and 4,1-NNL gave simple and very useful information for sequencing peptides and proteins in the MALDI-ISD experiments [2]. The observable fragment ions such as $[c+2H]^+$, $[z+2H]^+$, $[y+2H]^+$, $[a]^+$, $[x]^+$ and $[y]^+$ ions were exactly controlled by the sites of radical (unpaired electron) and Arg residue, and the RDD reactions obeyed simple and successive radical-initiative process.

Novel Aspect

It was indicated that $[a]^+$ and $[y_m+2H]^+$ ions were originated from beta-carbon-centered radical ions $[M-H+H]^{\cdot}$ and radical intermediates $[z_{m+1}+H]^{\cdot}$, respectively.

References

1. M.Takayama *et al.*, *Int. J. Mol. Sci.*, 15, 8428-8442 (2014).
2. K.Nagoshiet *al.*, *J. Am. Soc. Mass Spectrom.*, in press.

TOr.08-5 JUGGLING FRAGILE VAN DER WAALS COMPLEXES IN AN ION TRAP FOR MULTIPLEXED INFRARED PHOTODISSOCIATION

Larry Tesler - Matthew Bell - Taylor Harmon - Laura Bailey - Nicolas Polfer

University of Florida, University of Florida/Chemistry, Gainesville

Keywords: infrared photodissociation, cryogenic trap, multiplexing, spectroscopy, metabolites

Introduction

Cryogenic infrared spectra can in principle serve as a gold standard for molecular identification in MS [1]. However, a key failing of action spectroscopy remains its extremely low duty cycle, where only one analyte is probed at one wavelength step at one time. Here, we discuss the first mass-selective cryogenic trap [2], which allows the IR spectra of multiple analytes to be probed in parallel, given the predictable mass increase/decrease of tagging spectroscopy [3].

Methods

A custom cryogenic 2D linear ion trap (cryoLIT) is composed of rectilinear electrodes equipped with slits in the left/right rf electrodes for radial ejection. Temperature control is enabled via a helium cycle cryostat, cartridge heaters, and temperature sensors. A gas pulse of He (and 5% N₂) allows ion cooling and tagging. Ions are mass manipulated using an auxiliary dipolar waveform on the left/right rf electrodes that is resonant with the secular frequency.

Results

When cooling the cryoLIT to cryogenic temperatures, ions exit the trap at lower rf amplitudes in the mass instability scan, which is due to a shrinking of the trap. The magnitude of this shift is related to the expansion coefficient of the material, 0.33% for 17-PH steel and 0.65% for OFHC copper.

Two types of tagged ions will be considered here: 1) solvent-tagged ions (H₂O and CH₃CN) from the electrospray source, or 2) van der Waals N₂-tagged ions generated in the cryogenic trap. Given their extremely low binding energies, the parameters for generating and manipulating N₂-tagged ions are critical. Still, with optimized trap temperature (24K), pumpdown delays, and isolation waveforms a background-free mass isolation of N₂-tagged ions is feasible. High-quality IR spectra of small analytes (tryptophan, MDMA, Leu-enkephalin, etc.) can be recorded in this way and are very comparable to literature cryogenic IR spectra [4]. Current efforts focus on neon tagging, which may be the optimal tag in terms of mass increase (+20) and minimal IR spectral distortion.

Conclusions

The mass-selective cryoLIT is capable of measuring high-resolution IR spectra of multiple analytes in a multiplexed fashion at reasonable speeds (< 30 min), and at relevant concentrations (nM). This sets the stage for applying cryogenic IR spectroscopy as a bioanalytical tool for small molecule identification (e.g. metabolites). Two-trap designs can further improve speed by decoupling the pumpdown delay (100 ms) from mass manipulation and analysis (100 ms).

Novel Aspect

First mass-selective cryogenic 2D ion trap where extremely fragile van der Waals-tagged ions can

be mass isolated for background-free infrared photodissociation.

References

1. Cismesia AP, Bailey LS, Bell MR, Tesler LF, Polfer NC, *J. Am. Soc. Mass Spectrom.*, 27, 757 (2016).
2. Cismesia AP, Tesler LF, Bell MR, Bailey LS, Polfer NC, *J. Mass Spectrom.*, 52, 720 (2017).
3. Kamrath MZ, Relph RA, Guasco TL, Leavitt CM, Johnson MA, *Int. J. Mass Spectrom.*, 300, 91 (2011).
4. Pereverzev AY, Cheng X, Nagarnova NS, Reese DL, Steele RP, Boyarkin OV, *J. Phys. Chem A*, 120, 5598 (2016).

TO.09 - Metal and non-metal clusters: gas phase structures and reactivity

Chairs: Evan Williams, Pietro Traldi

TO.09-1 **Keynote: TANDEM CRYO ION TRAP MASS SPECTROMETRY FOR KINETICS AND SPECTROSCOPY OF TRANSITION METAL CLUSTERS**

Gereon Niedner-Schatteburg

Fachbereich Chemie and Forschungszentrum OPTIMAS, TU Kaiserslautern, Kaiserslautern, Germany

We have implemented cryo ion trap FT-ICR instrumentation with temperature control to as low as 11 Kelvin. With two such instruments we investigate the magnetism and the adsorption kinetics of transition metal clusters and complexes. We have upgraded one of these instruments by a cryo RF trap for tandem cryo operation. Spectroscopy comes into play through application of Infrared Photon Dissociation (IR-PD) by optical parametric oscillator/amplifier (OPO/OPA) photon sources. This enables one and two colour investigations of metal organic complexes^[1] and of transition metal cluster adsorbate complexes^[2-3].

Stepwise cryo adsorption kinetics of N₂ on Ni_n⁺ clusters reveal clearly discernible mono layer like adsorbate shells which interpret in the light of IR-PD spectra of the product complexes Ni_n(N₂)_m⁺. It emerges a characterization of the cluster surface morphology and its adsorbate induced reorganization^[4-5]. The according data on Fe_n⁺ cluster N₂ adsorption are subject of current interpretation.

Case studies on various metallophthalocyanine complexes found interpretation in terms of adsorbate induced spin quenching^[6], and recent findings of non classical blue shifts of N₂ stretching modes of adsorbates to Fe₃O(OAc)₆⁺ complexes received interpretation in terms of orbital interactions^[7].

The tandem cryo trap option comes to live by studying consecutive co-adsorptions of N₂ and H₂ to selected Ru₈⁺ clusters – as compared to the reverse order co-adsorption. Our findings allowed to elucidate the mutual interdependencies amongst the adsorption kinetics, their activation and migration^[8].

This work was supported by the Deutsche Forschungsgemeinschaft DFG in the framework of the transregional collaborative research center SFB/TRR 88 3MET.de.

References

1. J. Lang, M. Gaffga, F. Menges, GNS, *Two-color delay dependent IR probing of torsional isomerization in a AgL₁L₂⁺ complex*, *PCCP* **2014**, 16, 17417-17421.
2. S. Dillinger, J. Mohrbach, J. Hewer, M. Gaffga, GNS, *IR spectroscopy of N₂ adsorption on size selected cobalt cluster cations in isolation*, *PCCP* **2015**, 17, 10358-10362.
3. J. Mohrbach, S. Dillinger, GNS, *Cryo Kinetics and Spectroscopy of Cationic Nickel Clusters: Rough and Smooth Surfaces*, *J. Phys. Chem. C* **2017**, 121, 10907-10918.

4. J. Mohrbach, S. Dillinger, and GNS, *Probing cluster surface morphology by cryo kinetics of N₂ on cationic nickel clusters*, *J. Chem. Phys.* **2017**, *147*, 184304.
5. S. Dillinger, J. Mohrbach, and GNS, *Probing cluster surface morphology by cryo spectroscopy of N₂ on cationic nickel clusters*, *J. Chem. Phys.* **2017**, *147*, 184305.
6. S. Dillinger, J. Lang, and GNS, *Cryo IR Spectroscopy of Hemin⁺ Complexes in Isolation*, *J. Phys. Chem. A* **2017**, *121*, 7191-7196.
7. J. Lang, J. Mohrbach, S. Dillinger, J. M. Hewer, and GNS, *Vibrational blue shift of coordinated N₂ in Fe₃O(OAc)₆(N₂)_n⁺: "non-classical" dinitrogen complexes*, *Chem. Comm.* **2017**, *53*, 420-423.
8. S. Dillinger, M. P. Klein, A. Steiner, D. C. McDonald, M. A. Duncan, M. M. Kappes, and GNS, *Cryo IR Spectroscopy of N₂ and H₂ on Ru₈⁺: The Effect of N₂ on the H-Migration*, *J. Phys. Chem. Lett.* **2018**, *9*, 914-918.

TOr.09-2 UV-VIS SPECTROSCOPY AND PHOTOCHEMISTRY OF TRAPPED IONS USING FT-ICR-MS

Christian van der Linde - Milan Ončák - Nina K. Bersenkovitsch - Andreas Herburger - Thomas Taxer - Martin K. Beyer

Universität Innsbruck, Institut für Ionenphysik und Angewandte Physik, Innsbruck

Keywords: spectroscopy, photochemistry, FT-ICR-MS, atmosphere, magnesium cluster

Introduction

While infrared spectroscopy of ions is a well-established field, the potential of UV-VIS spectroscopy of ions is barely touched. The field was pioneered by Ben Freiser in the 1980s, using lamps and monochromators as light source.[1] Nowadays tunable optical parametric oscillators cover a wide range of wavelengths from UV to infrared which can be used for spectroscopy, e.g. the photodissociation of peptides.[2]

Methods

Two FT-ICR-MS setups are used: A modified CMSX47 instrument equipped with a 4.7T magnet, APEX II data station and a laser vaporization ion source and a Bruker Apex Qe, equipped with a nanobay console, electrospray ion source and a 9.4T magnet. Light from tunable OPO (EKSPLA NT342B) is guided into the FT-ICR cell for studying spectroscopy and photochemistry. Experiments are supported by theoretical calculations.

Results

FT-ICR-MS is, compared to time of flight techniques or quadrupole mass filters, inherently slow, with a duty cycle in the range of seconds. High-resolution spectroscopy over a wide wavelength range, which requires the recording of 10³-10⁴ mass spectra, is not the prime application of FT-MS instruments. However, FT-MS offers a combination of high mass resolution, long trapping times, flexible mass selection capabilities in a collision-free environment and simultaneous detection of all fragments with high fidelity of absolute intensities. These features gives FT-MS significant advantages over other techniques for specific problems.

Selected examples which illustrate these advantages include photochemical hydrogen formation in the Mg⁺(H₂O)_n, n = 1-5, system, or the spectroscopy of the atmospherically relevant CO₃⁻(H₂O)_{1,2} ion in combination with a nitrogen cooled FT-ICR cell.

Conclusions

The results show that FT-ICR-MS is a powerful tool in combination with theory for investigation of the photochemistry of cluster systems.

Novel Aspect

Interplay between single and two photon processes play an important role in photodissociation

of hydrated positively charged magnesium clusters.

References

1. Cassady C. J., Freiser B. S. *Journal of the American Chemical Society*, 106, 6176 (1984).
2. Herburger A., van der Linde C., Beyer M. K. *Physical Chemistry Chemical Physics*, 19, 10786 (2017).

TO.r.09-3 **ESI IN THE STUDY OF THE INTERACTIONS BETWEEN PHOSPHINO-CU(I) COMPLEXES AND CU(II) SALTS WITH MODEL PEPTIDES RELEVANT TO HUMAN COPPER TRANSPORTERS (HCTR)**

Francesco Tisato⁽¹⁾ - **Karolina Bossak-Ahmed**⁽²⁾ - **Marina Porchia**⁽¹⁾ - **Wojciech Bal**⁽²⁾ - **Arkadiusz Bonna**⁽³⁾

⁽¹⁾ ICMATE – CNR, Corso Stati Uniti, 4 – 35127 Padova, Italy - ⁽²⁾ Institute of Biochemistry and Biophysics, Polish Academy of Sciences, Pawińskiego 5a, 02-106 Warsaw, Poland - ⁽³⁾ Department of Biochemistry, University of Cambridge, Tennis Court Road, Downing Site, Cambridge CB2 1QW, United Kingdom

Keywords: Copper anticancer drugs, internalization of copper in cancer cells, Cu(I) and Cu(II), model tripeptides, human Copper Transporter1 (hCtr1)

Introduction

Human Copper Transporter 1 (hCtr1) is a 190 amino-acids trans-membrane protein mainly devoted to the recruitment of metal ions from the extra-cellular fluids [1]. The extracellular domain of hCtr1 contains histidine-rich and methionine-rich sequences likely responsible of metal ions recruitment. The terminal portion of the protein extruding from the cell membrane contains the characteristic Met-Asp-His peptide sequence.

Methods

The interaction of Cu(I), Ag(I) and Cu(II) ions comprised in appropriate precursors (*i.e.* $[\text{Cu}(\text{P})_4]^+$, $[\text{Ag}(\text{P})_4]^+$ and CuCl_2 (P = tertiary phosphine)) with selected model tripeptides XYHis (X = Gly, Met; Y = Gly, Asn, Asp) and some polypeptides relevant to hCtr1 has been investigated by electrospray ionization mass spectrometry, and the interaction products have been characterized by multiple collisional experiments, using an ion trap mass instrument.

Results

Model XYHis tripeptides starting from the designed GlyGlyHis derivative of the ATCUN family [2] are able to mine Cu(I) and Ag(I) from $[\text{M}(\text{I})(\text{P})_4]^+$ inorganic precursors through the formation of the $[\text{M}(\text{I})(\text{XYHis})]^+$ adducts under ESI(+)-MS conditions [3]. Also Cu(II) is wrapped by XYHis tripeptides via $[\text{Cu}(\text{II})(\text{XYHis}) - \text{H}]^+$ adducts. MSⁿ of these metal-tripeptide complexes indicate that His is bound strongly to the metal in all cases. While CID pathways of $[\text{M}(\text{I})(\text{XYHis})]^+$ ions proceeded through losses of neutral fragments, MSⁿ fragmentation of $[\text{Cu}(\text{II})(\text{XYHis}) - \text{H}]^+$ ions instead took place mainly through intramolecular electron transfer reactions including the reduction of Cu(II) to Cu(I) and the formation of fragment radical cations. The interaction of $[\text{M}(\text{I})(\text{P})_4]^+$ precursors with polypeptide sequences representing the terminal extra-cellular portion of hCtr1 (MetAspHis, MetAspHisSerHisHis and MetAspHisSerHisHisMetGlyMetSer) confirms the ability of these peptides to mine and recruit M(I) ions from inorganic precursors.

Conclusions

Model tripeptides having a free-amino group at N-terminus and either a free carboxylic or an amidated group at the C-terminus of His in third position, and polypeptide sequences relevant to hCtr1 are able to mine Cu(I), Ag(I) and Cu(II) from inorganic precursors via formation of $[\text{M}(\text{tri/polypeptide})]^+$ adducts. This experimental evidence is significant to shed light in the cell internalization process of anticancer phosphino-Cu(I) drugs.

Novel Aspect

The observation that Cu(II) comprised in the $[\text{Cu(II)(XYHis)} - \text{H}]^+$ ions undergoes redox reaction under ESI(+)MSⁿ conditions might be relevant in the claimed reduction processes of physiological Cu(II).

References

1. Lee J., Thiele D.J., et al., *J. Biol. Chem.*, 277, 4380 (2002).
2. Harford C., Sarkar, B., *Acc. Chem. Res.* 30, 123 (1997)
3. Tisato F., Traldi, P., et al. *Eur. J. Mass Spectro.* 22, 275 (2016)

TOr.09-4 STRUCTURE AND REACTIVITY OF ANTIMONATO POLYOXOVANATES IN SOLUTION AND IN THE GAS PHASE

Ulrike Warzok⁽¹⁾ - Lisa K. Mahnke⁽²⁾ - Michael Wendt⁽²⁾ - Wolfgang Bensch⁽²⁾ - Christoph A. Schalley⁽¹⁾

⁽¹⁾ Institute for Chemistry and Biochemistry, Free University Berlin, Berlin - ⁽²⁾ Institute for Inorganic Chemistry, Christian-Albrechts-Universität zu Kiel, Kiel

Keywords: Antimonato polyoxovanadates, solution and gas phase reactivity, isotope exchange experiments, CID, mechanistic studies

Introduction

Polyoxometalates cover a broad spectrum of interesting properties related to the nature and number of incorporated metals, their size and shape, their redox behavior, magnetic and photochemical characteristics. Possible applications are manifold. Yet to harness a certain functionality, comprehensive insights into their structure and reactivity are crucial.

Methods

We are presenting the investigation of structure and reactivity in solution and in the gas phase of different antimonato polyoxovanadates by electrospray ionization mass spectrometry (ESI-MS), isotope exchange experiments in solution, collision-induced dissociation tandem MS and ion mobility MS experiments conducted at a SYNAPT G2-Si (Waters, Manchester).

Results

The first water-soluble antimonato polyoxovanadate cluster was demonstrated to be present in solution as two distinct species, which both possess an intact, closed cage-like $\{\text{V}_{15}\text{Sb}_6\}$ structure: a water-encapsulating and a water-free cluster core. Its chemical behavior is greatly dependent on the encapsulated water molecule, since a transduction of inner-phase reactivity of the guest results in changes in the outer-phase reactivity of the cluster cage.[1] Moreover, the surrounding cationic metal-complexes and the ligand environment, as well as additives influence the reactivity of the $\{\text{V}_{15}\text{Sb}_6\}$ cluster in solution.[2] In this regard, a $\{\text{Sb}_6\text{V}_{15}\} \rightarrow \{\text{Sb}_8\text{V}_{14}\}$ cluster transition was discovered by ESI-MS. This reaction enables in combination with a change in the ligand sphere the formation of an unusual, meta-stable $\alpha_1^* - \{\text{V}_{14}\text{Sb}_8\}$ isomer which forms dimeric super structures. These were investigated in the solid state, in solution and in the gas phase.[3]

Conclusions

ESI-MS and IMS could serve as a probe for the elucidation of solution and gas-phase reactivity of metal-oxo clusters, as well as the identification and structural analysis of new products. These methods can provide the rational guideline for the synthesis of new compounds and for the assembly of resulting complex supramolecular architectures.

Novel Aspect

We present the first comprehensive study on the structure as well as solution and gas-phase reactivity of antimonato polyoxovanadate clusters.

References

1. M. Wendt, U. Warzok, C. Näther, J. van Leusen, P. Kögerler, C. A. Schalley, W. Bensch, *Chem. Sci.* **2016**, 7, 2684–2694.
2. Mahnke, U. Warzok, M. Lin, C. Näther, C.A. Schalley, W. Bensch, *Chem. Eur. J.* **2018**, 24, 5522–5528.
3. L. K. Mahnke, A. Kondinski, U. Warzok, C. Näther, J. van Leusen, C.A. Schalley, K.Y. Monakhov, P. Kögerler, W. Bensch, *Angew. Chem. Int. Ed.* **2018**, 57, 2972–2975.

TO.r.09-5 FT-ICR MS STUDIES OF THIOMOLYBDATE CLUSTERS AS MODEL CATALYSTS FOR HYDROGEN EVOLUTION

Aristeidis Baloglou⁽¹⁾ - **Milan Ončák**⁽¹⁾ - **Christian van der Linde**⁽¹⁾ - **Marie-Luise Grutza**⁽²⁾ - **Philipp Kurz**⁽²⁾ - **Martin K. Beyer**⁽¹⁾

⁽¹⁾ Universität Innsbruck, Institut für Ionenphysik und Angewandte Physik, Innsbruck - ⁽²⁾ Albert-Ludwigs-Universität Freiburg, Institut für Anorganische und Analytische Chemie, Freiburg

Keywords: molybdenum sulfide; catalysis; hydrogen evolution; action spectroscopy; dissociation

Introduction

Molybdenum sulfide-based catalysts have the potential to replace platinum as a hydrogen evolution reaction (HER) catalyst in electrochemical water splitting [1]. A better understanding of the reaction mechanisms contributes to the targeted optimization of the catalysts. In this contribution, individual reaction steps are investigated quantitatively employing precisely defined gas-phase models, like $\text{Mo}_3\text{S}_{13}^{(2-)}$ and $\text{Mo}_2\text{S}_{12}^{(2-)}$.

Methods

Fourier transform ion cyclotron mass spectrometry (FT-ICR MS) is used, as it is an excellent tool to examine ion-molecule reactions in the gas phase [2], including catalytic cycles. By combining FT-ICR MS with IR and UV/VIS laser spectroscopy, the reactants, products and reactive intermediates can be further characterized.

Results

IR-spectroscopy experiments with the protonated species $\text{HMo}_3\text{S}_{13}^-$ show a distinct H-S stretching vibration. Comparing this result with density functional theory (DFT) calculations suggests that the proton is bound to an edge disulfide unit of the cluster. In gas phase, the cluster $\text{Mo}_3\text{S}_{13}^{(2-)}$ unfolds from the highly symmetric Mo_3 -ring structure, found in the solid state, to form irregular structures with a variety of mono- and disulfide ligands. Structural changes are also expected for molybdenum sulfide clusters deposited on electrodes. Investigation of the structural flexibility and chemical properties of the different sulfur moieties contributes to an understanding of the catalytic activity of molybdenum sulfide surfaces.

Conclusions

Structural features evolve in the gas phase that elucidate the conformational flexibility and relative stability of reactive molybdenum sulfide species. Hydrogen is bound at sulfur centers, as revealed by the infrared signature of the S-H bond. Electronic excitation spectra indicate the availability of low-lying excited states for some species generated by mild collisional activation.

Novel Aspect

Molybdenum sulfide clusters are synthesized with isotopically enriched molybdenum and characterized by mass spectrometry in combination with optical excitation.

References

- [1] Karunadasa H. I., Montalvo E., Sun Y., Majda M., Long J. R., Chang C. J., *Science*, 335, 698 (2012).
- [2] Baloglou A., Ončák M., van der Linde C., Beyer M. K., *Top. Catal.* 61, 20 (2017).

TOr.10 – Astrochemistry

Chairs: Paul Mayer, Simonetta Fornarini

TOr.10-1 **Keynote: COSMIC PAHS: FORMATION AND UV PROCESSING**

Christine Joblin

Institut de Recherche en Astrophysique et Planétologie (IRAP), Université de Toulouse, CNRS, CNES, Toulouse

Keywords: Laboratory astrophysics, polycyclic aromatic hydrocarbons, ion traps, UV processing, gas-phase formation

Introduction

Polycyclic aromatic hydrocarbons (PAHs) play a key role in the physical and chemical evolution of UV-irradiated astrophysical environments from protoplanetary disks to galaxies. In order to progress on our knowledge of these species, laboratory astrophysics activities have both to address the interaction of PAHs with UV photons, which is key to their stability, and investigate possible formation routes for these large molecules in cosmic environments.

Methods

The study of the photophysics of isolated PAHs, including ionisation, dissociation and radiative cooling benefits from different setups (molecular jets, ion traps, storage rings) and their coupling with VUV synchrotron radiation [1-5]. PAH formation mechanisms are investigated using gas-phase chemical reactors, in order to identify key reactions [6] or evaluate the impact of various experimental parameters on the production of carbon dust and PAHs [7-10].

Results

In the recent years, we have been able to better quantify the photophysics of isolated PAHs, including the dynamics of ionisation [11], as well as the rates of radiative cooling [5] and fragmentation [1]. The contribution of PAHs to the formation of H₂ has gained further support [12-14]. The detection of the fullerene C₆₀ in cosmic environments has led to considerations on the formation of this species by VUV photoprocessing of large PAHs [15].

On the other hand, activities are on going within the Nanocosmos ERC Synergy project [16] to study the formation of carbon dust in conditions that prevail in evolved stars, where this dust is likely formed. More specifically, the Stardust machine [10] aims at mimicking the key steps involved in dust formation around evolved stars. The AROMA setup was designed to support these experiments by analysing the molecular content of the synthesised dust [17]. This analysis, combined with that of meteorite samples, is expected to provide new insights into chemical pathways leading to the formation of cosmic PAHs.

Conclusions

Based on their photophysical properties, only large PAHs containing ~100 atoms are predicted to survive in the harsh conditions of UV-irradiated regions. On the other hand, these large PAHs do not appear to be formed in gas-phase reactors neither be present in meteorite samples.

Understanding this point requires a multidisciplinary approach at the frontier between astronomy, molecular physics, and chemistry. Mass spectrometry remains a key tool in these studies.

Novel Aspect

In order to identify relevant cosmic PAHs, dedicated laboratory setups are now available to study their photophysics and tackle the chemistry involved in their formation.

References

1. West, B., Rodriguez Castillo, S., Sit, A., Mohamad, S., Lowe B., Joblin, C., Bodi, A., Mayer, P. M., *Phys. Chem. Chem. Phys.*, 20, 7195 (2018)
2. Rouillé, G., Krasnokutski, S. A., Fulvio, D., Jäger, C., Henning, Th., Garcia, G. A., Tang, X.-F., Nahon, L., *Astrophys. J.*, 810, id. 114 (2015)
3. Zhen, J., Rodriguez Castillo, S., Joblin, C., Mulas, G., Sabbah, H., Giuliani, A., Nahon, L., Martin, S., Champeaux, J.-P., Mayer, P., *Astrophys. J.*, 822, id. 113 (2016)
4. Zhen, J., Castellanos, P., Paardekooper, D. M., Ligterink, N., Linnartz, H., Nahon, L., Joblin, C., Tielens, A.G.G.M., *Astrophys. J. Lett.*, 804, L7 (2015)
5. Martin, S., Ji, M., Bernard, J., Brédy, R., Concina, B., Allouche, A. R., Joblin, C., Ortega, C., Cassimi, A., Ngono-ravache, Y., and Chen L., *Phys. Rev. A*, 92, id.053425 (2015)
6. Zhao, L., Kaiser, R. I., Xu, B., Ablikim, U., Ahmed, M., Joshi, D., Veber, G., Fischer, F. R., Mebel, A. M., *Nature Astron.*, 2, 413 (2018)
7. Jäger, C., Huisken, F., Mutschke, H., Jansa, I. Llamas, Henning, Th., *Astrophys. J.*, 696, 706 (2009)
8. Biennier, L., Georges, R., Chandrasekaran, V., Rowe, B., Bataille, T., Jayaram, V., Reddy, K.P.J., Arunan, E., *Carbon*, 47,3295 (2009)
9. Contreras, C. S., Salama, F., *Astrophys. J. Supp.*, 208, id. 6 (2013)
10. Martínez, L., Lauwaet, K., Santoro, G., Sobrado, J. M., Peláez, R. J., Herrero, V. J., Tanarro, I., Ellis, G., Cernicharo, J., Joblin, C., Huttel, Y., Martín-Gago, J. Á., *Scientific Reports*, (2018)
11. Marciniak, A., Despré, V., Barillot, T., Rouzée, A., Galbraith, M.C.E., Klei, J., Yang, C.-H., Smeenk, C.T.L., Lorient, V., Reddy, S. Nagaprasad, Tielens, A.G.G.M., Mahapatra, S., Kuleff, A. I., Vrakking, M.J.J., Lépine, F., *Nature Comm.*, 6, 7909 (2015)
12. Paris, C., Alcamí, M., Martín, F., and Díaz-Tendero, S., *J. Chem. Phys.*, 140, 204307 (2014)
13. Chen, T., Gatchell, M., Stockett, M. H., Delaunay, R., Domaracka, A., Micelotta, E. R., Tielens, A. G. G. M., Rousseau, P., Adoui, L., Huber, B. A., Schmidt, H. T., Cederquist, H., Zettergren, H., *J. Chem. Phys.*, 142, 144305 (2015)
14. Rodriguez Castillo, S., Simon, A., Joblin, C., *Int. J. Mass Spectrom.*, doi:10.1016/j.ijms.2017.09.013 (2018)
15. Berné O., Montillaud J., Joblin C., *Astron. & Astrophys.*, 577, id.A133 (2015); Corrigendum *Astron. & Astrophys.*, 588, id.C1 (2016)
16. <https://nanocosmos.iff.csic.es/>
17. Sabbah, H., Bonnamy, A., Papanastasiou, D., Cernicharo, J., Martín-Gago, J.-A., Joblin, C., *Astrophys. J.*, 843, id. 34 (2017)

TOr.10-2 RELAXATION AND FRAGMENTATION DYNAMICS OF POLYCYCLIC AROMATIC HYDROCARBONS EXCITED BY AN ULTRASHORT XUV PULSE

Marius Herve⁽¹⁾ - **Pablo Castellanos**⁽²⁾ - **Gabriel Karras**⁽¹⁾ - **Victor Despre**⁽¹⁾ - **Alexandre Marciniak**⁽¹⁾ - **Eric Constant**⁽¹⁾ - **Vincent Lorient**⁽¹⁾ - **Alexander Kuleff**⁽³⁾ - **Alexander Tielens**⁽²⁾ - **Franck Lepine**⁽¹⁾

⁽¹⁾ Institut Lumière Matière (ILM), Université Lyon 1, Lyon - ⁽²⁾ Leiden Observatory, Leiden University, Leiden - ⁽³⁾ Physikalisch-Chemisches Institut, Universität Heidelberg, Heidelberg

Keywords: PAHs; XUV excitation; femto-astrochemistry; cation fragmentation dynamics; many-body quantum effects

Introduction

Sources of ultrafast XUV pulses have brought new possibilities of exploring relaxation dynamics of gas phase molecules following high-energy absorption. We used XUV excitation to produce highly excited states in Polycyclic Aromatic Hydrocarbons (PAHs) cations, thought as important components of interstellar media, and probed their subsequent femtosecond dynamics. This reveals that many-body quantum effects could have consequences in astrochemistry models.

Methods

We used an XUV pump-IR probe scheme coupled to a time-of-flight setup [1]. Femtosecond XUV pump pulses centered around 23 eV were generated by High-order Harmonics Generation (HHG), and recombined with delayed 800 nm probe pulses to interact with gas phase neutral PAHs of various size (from Naphtalene $C_{10}H_8$ to Hexabenzocoronene $C_{42}H_{18}$), creating XUV excited cations probed by IR photons. Fragments and dication yields were then measured as a function of the XUV-IR delay.

Results

In all the studied PAHs, dication yield reveals an instantaneous population of cationic states, followed by an ultrafast decay, ranging from 29 fs for naphtalene to 59 fs for hexabenzocoronene. XUV-induced shake-up ionization indeed forms highly correlated cationic states lying just below the second ionization potential, whose relaxation dynamics to lower states is governed by coupled electron and nuclear motion [2]. Addition of a delayed IR photon thus probes this dynamics by second ionization. Increase of this timescale with PAH size is explained by the higher density of shake-up states through which initial population has to decay. We further observed population of fragmentation channels on the 100 fs timescale (e.g. 119 fs for H loss of Pyrene) in the fragment yields. This longer timescale corresponds to the vibrational energy redistribution in the molecule following the non-adiabatic shake-up states relaxation, leading finally to fragmentation. It reveals the energy spread to nuclear degrees of freedom.

Conclusions

We reported dynamics occurring in small and large PAHs after XUV absorption. Relaxation dynamics of highly excited cations occurs on a similar timescale for all the PAHs, increasing with PAHs size, followed by energy redistribution to vibrations. It demonstrates that relaxation of these highly correlated states through many-body quantum effects is a general mechanism for PAHs. These general trends could be incorporated in the description of interstellar media [3].

Novel Aspect

Femtosecond XUV sources were used to measure ultrafast relaxation dynamics of astro-relevant PAHs, where electron and nuclei are fully coupled.

References

1. Lorient V. *et al*, J. Phys.: Conf. Ser., 635, 012006 (2015)
2. Marciniak A. *et al*, Nat. Commun., 6, 7909 (2015)
3. Tielens A., Ann. Rev. Astron. Astrophys., 46,289 (2008)

TOOr.10-3 EXPERIMENTAL STUDY OF HIGH-ENERGY ION/MOLECULE INTERACTIONS RELEVANT TO PLANETARY ATMOSPHERES

Vera Krizova ⁽¹⁾ - Jan Zabka ⁽²⁾ - Illia Zymak ⁽²⁾ - Miroslav Polasek ⁽²⁾

⁽¹⁾ Charles University in Prague, Faculty of Science, Prague - ⁽²⁾ J. Heyrovsky Institute of Physical Chemistry of the CAS, Department of Chemistry of Ions in Gaseous Phase, Prague

Keywords: ion/molecule interactions, high-energy interactions, ion-induced processes, planetary atmospheres

Introduction

Atmospheres of Titan and Jupiter have already been investigated by several probes in the past. Atmosphere of Titan is composed of nitrogen and methane and other minor constituents [1], while Jupiter's atmosphere comprises mainly hydrogen and helium in addition to other less abundant compounds [2]. Constituents of both atmospheres undergo ionization by energetic heavy ions present in magnetospheres of Saturn [3] and Jupiter [4] respectively.

Methods

The experiments were performed on a modified sector-type mass spectrometer having a quadrupole analyzer attached to a collision chamber located between magnetic and electrostatic sector. The incident ions were formed in EI source, accelerated to 5 – 10 keV, mass-selected by the magnetic sector and directed through the collision chamber. The secondary ions stemming from the neutrals admitted to the collision chamber were analyzed by the quadrupole mass analyzer.

Results

In this study we focused on projectiles and neutral targets relevant to ionospheres of Titan and Jupiter. When projectiles interact with a neutral target, two processes can occur: ionization and electron transfer. Occurrence of one or the other depends on thermochemistry of the reaction. Objective of this study is to investigate the amount of internal energy deposited in nascent secondary ions which is reflected by the extent of fragmentation of the ionized neutral. It was found out that in the case of atomic projectiles, fragmentation of the ionized neutral strongly depends on the ionization energy of the projectile but depends very little on the center-of-mass collision energy in the collision energy range used. However in the case of molecular projectiles, the extent of fragmentation is generally very low, i.e. only molecular ions are formed, and no apparent dependencies of ionization energies, sizes or chemical structures of ionic projectiles were found.

Conclusions

For atomic projectiles, it was found that the extent of fragmentation depends on thermochemistry of electron transfer: the higher difference in ionization energies, the higher the extent of fragmentation is. For molecular projectiles the extent of fragmentation seems to be inversely proportional to the size of the projectile.

Novel Aspect

To our knowledge, processes investigated in this study with the particular sets of projectiles and neutrals have not yet been investigated experimentally.

Acknowledgment

The study was supported by the Charles University, project GA UK No. 928417.

References

1. Coustenis A., *Space Sci. Rev.* **116**, 171 (2005)
2. Niemann H. B., et al., *J. Geophys. Res.* **103**, 22831 (1998)
3. Arridge C. S., et al., *Space Sci. Rev.* **162**, 25 (2011)
4. Frank L. A., et al., *Space Sci. Rev.* **60**, 283 (1992)

TOr.10-4 INVESTIGATING THE POSSIBILITY OF FORMING PREBIOTIC MOLECULES BY GAS PHASE ION-MOLECULE REACTIONS. APPLICATION TO GLYCINE

Riccardo Spezia⁽¹⁾ - ***Yannick Jeanvoine***⁽²⁾ - ***Antonio Largo***⁽³⁾ - ***William L. Hase***⁽⁴⁾ - ***Joel Lemaire***⁽⁵⁾ - ***Debora Scuderi***⁽⁶⁾

⁽¹⁾ CNRS, Sorbonne Université, Laboratoire de Chimie Théorique, Paris - ⁽²⁾ Université d'Evry, Laboratoire Analyse et Modélisation pour la Biologie et l'Environnement, Evry - ⁽³⁾ Universidad de Valladolid, Departamento de Química Física, Valladolid - ⁽⁴⁾ Texas Tech University, Department of Chemistry and Biochemistry, Lubbock - ⁽⁵⁾ CNRS, Université Paris Sud, Laboratoire de Chimie Physique, Orsay - ⁽⁶⁾ Université Paris Sud, Laboratoire de Chimie Physique, Orsay

Keywords: Ion-molecule reactions; prebiotic chemistry; chemical dynamics simulations; Infra-red multi photon dissociation (IRMPD)

Introduction

In last years, progresses in astrophysics have revealed the presence in the interstellar medium (ISM) and in comets of many organic molecules, like glycine [1], urea [2] and formamide [3]. Due to the extreme conditions in space, it is not clear how such prebiotic molecules can be synthesized. Gas-phase ion chemistry, combining simulations and experiments, can propose some synthetic routes.

Methods

Chemical dynamics simulations and quantum chemistry were used to theoretically investigate ion-molecule reactions [4]. Simulations were done at different collision energies. Experimentally, we have used a Paul ion trap where we have introduced the ions and the neutrals. IR laser was used to characterize the products [5] and to induce the reactivity of the stabilized complexes [6].

Results

To form glycine, we have considered the reaction: $\text{NH}_3\text{OH}^+ + \text{CH}_3\text{COOH}$. Bohme and co-workers reported that it produces protonated glycine in a flow-tube [7].

Chemical dynamics simulations, tuning collision energy, form an ion with m/z 76 and water. However, the ions formed do not have the structure of protonated glycine [8].

Experimentally, the same reaction done in the ion trap conditions does not form the ion m/z 76, but the ion m/z 94, corresponding to the complex made by NH_3OH^+ and CH_3COOH . It was irradiated by IR and characterized comparing IRMPD and theoretical spectra.

When doing IRMPD on the m/z 94, one of the photo-fragment is the ion m/z 76, corresponding to protonated glycine, suggesting a possible role of IR radiation in the synthesis of such prebiotic molecule.

Preliminary calculations [9] and experiments on urea formation suggests that a similar picture can occur also in this case.

Conclusions

We have shown that glycine can be formed by ion molecule collisions between protonated hydroxylamine and methanol. This reaction can occur if some translational energy is given or if the formed complex is irradiated by an IR laser.

The species observed in simulations have the same m/z of protonated glycine but different structures, suggesting that radio-astronomy should consider this possibility.

Novel Aspect

Simulations and experiments of ion-molecule reactions suggest new synthetic route to form prebiotic molecules under astrophysical conditions.

References

1. Altwegg, K., et al., *Sci. Adv.*, 2, e1600385 (2016).
2. Remijan, A. J., Snyder, L. E., McGuire, B. A., et al., *Astrophys. J.*, 783, 77 (2014).

3. Rubin, R. H., Swenson, G. W., Jr., Solomon, R. C., Flygare, H. L., *Astrophys. J. Lett.*, 169, L39 (1971).
4. Spezia, R., Jeanvoine, Y., Hase, W.L., Song, K., Largo, A., *Astrophys. J.*, 826, 107 (2016).
5. Mac Aleese, L., Simon, A., McMahon, T.B., Ortega, J.M., Scuderi, D., Lemaire, J., Maitre, P., *Int. J. Mass Spectrom.* 249, 14 (2006)
6. Scuderi, D., Lepere, V., Piani, G., Bouchet, A., Zehnacker-Rentien, A., *J. Phys. Chem. Lett.* 5, 56 (2014).
7. Snow, L.J., Orlova, G., Blagojevic, V., Bohme, D. K., *J. Am. Chem. Soc.*, 129, 9910 (2007).
8. Jeanvoine, Y., Largo, A., Hase, W.L., Spezia, R., *J. Phys. Chem. A*, 122, 869 (2018)
9. Siro Brigiano, F., Jeanvoine, Y., Largo A., Spezia, R., *Astron. Astrophys.*, in press. DOI: 10.1051/0004-6361/201731610

TOr.10-5 THE THRESHOLD DISPLACEMENT ENERGY OF BUCKMINSTERFULLERENE C₆₀ AND FORMATION OF THE ENDOHEDRAL DEFECT FULLERENE HE@C₅₉

Mark Stockett⁽¹⁾ - Michael Wolf⁽¹⁾ - Michael Gatchell⁽²⁾ - Henning Schmidt⁽¹⁾ - Henning Zettergren⁽¹⁾ - Henrik Cederquist⁽¹⁾

⁽¹⁾ Stockholm University, Department of Physics, Stockholm - ⁽²⁾ Universität Innsbruck, Institut für Ionenphysik und Angewandte Physik, Innsbruck

Keywords: fullerenes, collisions, knockout, endohedrals, supernovae

Introduction

The observation of fullerenes in various astronomical environments, and the identification of C₆₀⁺ as the carrier of several Diffuse Interstellar Bands, has lead to a surge in interest in the formation and destruction mechanisms of fullerenes under different astrophysical conditions. Non-statistical fragmentation of fullerenes and other large molecules following collisions with energetic H and He atoms in supernova shocks is considered a key mechanism [1].

Methods

C₆₀⁻ ions are produced by electrospray, selected with a quadrupole mass filter, and accelerated to 1-15 keV. Collisions with He at center-of-mass energies of 20-80 eV mimic supernova processing. Negatively and positively charged product distributions are measured with an electrostatic energy analyzer. Classical molecular dynamics (MD) simulations enable detailed analysis of collisions.

Results

As the electron affinity of C₆₀ is much lower than any of its dissociation energies, collisions depositing enough energy to induce statistical unimolecular dissociation most likely lead to electron loss and thus do not contribute to the negative ion product spectrum. Non-statistical fragmentation *i.e.* single carbon knockout is observed directly by detecting C₅₉⁻. We observe for the first time the endohedral defect fullerene complex He@C₅₉⁻. Secondary decay products of these ions, C₅₈⁻ and He@C₅₈⁻, are also observed. We measure the threshold center-of-mass energy for knocking out a single carbon atom from C₆₀ in collisions with He. Combining this with our MD simulations, we determine a semi-empirical value for the (projectile-independent) threshold displacement energy, the minimum energy needed to remove a single carbon atom from the C₆₀ cage [2]. Our value, 24.1(5) eV, is much higher than generally assumed previously for fullerenes (around 15 eV), and is similar to that determined for graphene [3] and Polycyclic Aromatic Hydrocarbons [4].

Conclusions

Non-statistical fragmentation of fullerenes and other large carbonaceous molecules is an

important destruction mechanism in collisions like those taking place in supernova shocks. The products of these collisions are surprisingly robust, and may serve as reactive seeds for bottom-up growth of larger molecules in the Interstellar Medium. The observation of $\text{He}@C_{59}^-$ is a testament to the complexity of fullerene reactions.

Novel Aspect

We report the first observation of endohedral defect fullerene ions and the first measurement of the threshold displacement energy of a fullerene.

References

1. Micelotta E. R. , Jones A. P., Tielens A. G. G. M., *Astron. Astrophys.* 510, A36 (2010).
2. Stockett M. H. *et al.*, <http://arxiv.org/abs/1804.04374>
3. Meyer J. C. *et al.*, *Phys. Rev. Lett.* 108, 196102 (2012).
4. Stockett M. H. *et al.*, *J. Phys. Chem. Lett.* 6, 4504 (2015).

WEDNESDAY, August 29th, 2018

PL04: Plenary Lecture: Curt Brunnée Medalist

Chair: Scott A. McLuckey

PL04 LITHOGRAPHICALLY PATTERNED ELECTRODES FOR MINIATURIZED ION TRAP MASS SPECTROMETERS AND OTHER ION OPTICS DEVICES

Daniel E. Austin - Stephen A. Lammert - Aaron R. Hawkins - Yuan Tian - Ailin Li - Trevor Decker - Brett Hansen - Qinghao Wu

Brigham Young University, Department of Chemistry and Biochemistry and Department of Electrical and Computer Engineering, Provo (USA)

Keywords: microfabrication, ion trap, planar electrode, charge detection

Introduction

All mass and mobility analyzers rely on carefully shaped electric fields (sometimes combined with magnetic fields) to separate ions. These electric fields are traditionally established using shaped metal electrodes. We have demonstrated several types of devices in which electrodes are patterned onto circuitboard or onto rigid substrates (ceramic, glass). Advantages include miniaturization, novel devices, and complex electrode arrangements.

Methods

We have developed quadrupole and linear ion traps, concentric dual ion traps, multi-stage image-charge detectors, charge detection mass analyzers, and electrostatic ion beam traps. Recent studies have focused on substrate misalignment effects, modeling the electric field, and microparticle analysis.

Results

We have demonstrated the operation of a sub-mm linear ion trap with better than unit mass resolution for small organic molecules. We have also observed the effect of misalignment between electrodes in several degrees of freedom (corresponding to the possible orientations of the substrates), and compared the measurements to simulations of misalignment. Using angled electrodes printed on circuit board material, we can measure the trajectory of microparticles

passing through a vacuum system, such as the divergence from an aerodynamic lens or microparticle-surface impact and rebound. We are also developing a microparticle analyzer (mass, charge) for a future mission to study aerosols on Mars.

Conclusions

These devices demonstrate the utility and flexibility of creating complex and miniaturized devices for ion analysis using electrode printing on substrates. Electrode alignment, which is critical for many applications, is easily explored using alignment between plates. Charge detection with printed electrodes has the same sensitivity as conventional metal electrodes.

Novel Aspect

Novel approach to making ion optics devices including ion traps and charge detection analyzers. First study of misalignment effects in ion traps in multiple degrees of freedom.

WOr.01 - Native MS and structural proteomics

Chairs: Frank Sobott, Nathalie Norais

WOr.01-1 **Keynote: NATIVE MASS SPECTROMETRY: A TECHNIQUE OF MANY FACES**

Michal Sharon - Gili Ben-Nissan - Shay Vimer - Jelena Cveticanin

Department of Biomolecular Sciences, Weizmann Institute of Science, Rehovot 7610001, Israel

Keywords: Native MS, Protein Complexes, Antibodies, Crude Proteins Samples, Protein-Ligand interactions

Introduction

Key processes in a cell are carried out by protein complexes whose functions require tight regulation. The classical, high-resolution methods used for understanding protein regulation are often insufficient when attempting to study large protein assemblies, due to their size, heterogeneity and dynamic composition – thus, new methodologies are essential. Native MS is one such method.

Methods

I will present our studies involving different developments and applications of native MS approaches with examples focusing on protein-protein interactions, protein-ligand systems and antibodies^{3,4}. I will also show a direct MS method that enables characterization of recombinant proteins directly from crude samples. Involving proteins produced in bacteria, yeast, insect and human cells.

Results

The direct-MS approach can be carried out for proteins produced in prokaryotic and eukaryotic expression systems. The method provides immediate information on the expression, identity, solubility, oligomeric state, post-translational modification, ligand binding, overall structure and stability of overproduced proteins. I will show that the expression profile, post-translational modifications, overall structure and stability of engineered therapeutic antibodies, as well as glycosylated proteoforms of a modified protein can be immediately assessed by analysis of the crude samples. Moreover, I will demonstrate that the method is applicable for throughput characterization of protein variants, as mutants and engineered constructs. Our method is general and can be applied on various mass spectrometry platforms.

Conclusions

The major advantage of the method is that it considerably shortens the time gap between protein production and characterization while saving the burden associated with purifying inactive variants. The approach is particularly suitable for selecting the ideal protein expression host system, optimizing codon usage, expression levels and quality and employing high throughput investigations of lead candidates.

Novel Aspect

Given the broad use of recombinant proteins in biochemical investigations and in industrial and therapeutic applications the method is expected to be widely utilized.

References

1. Sharon, M. *Science*. 340,1059-1060, (2013).
2. Dyachenko, A., Gruber R., Shimon, L., Horovitz, A. Sharon, M., *PNAS*. 110, 7235-7239 (2013).
3. Gan, J., Ben-Nissan, G., Arkind, G., Tarnavsky, M., Trudeau, D., Noda Garcia, L., Tawfik, D., Sharon, M., *Anal. Chem.* 89, 4398-4404,(2017).
4. Ben-Nissan, G., Belov, M. E., Morgenstern, D., Levin, Y., Dym, O., Arkind, G., Lipson, C., Makarov, A. A., Sharon, M.,*Anal. Chem.*89, 4708-4715, (2017).

WOr.01-2 NATIVE MASS SPECTROMETRY REVEAL FUNCTIONING MECHANISM OF MEMBRANE PROTEIN TRANSLOCATION PATHWAY

Shahid Mehmood⁽¹⁾ - **Melanie McDowell**⁽²⁾ - **Irmgard Sinning**⁽²⁾ - **Carol V. Robinson**⁽¹⁾

⁽¹⁾ Department of Chemistry, University of Oxford, Oxford - ⁽²⁾ Biochemie Zentrum, Heidelberg University, Heidelberg

Keywords: Membrane proteins; Tail-anchored protein translocation, native mass spectrometry;

Introduction

Native mass spectrometry (MS) has a strong potential for its application towards the systems where other biophysical techniques are futile (1). In yeast, the dedicated guided entry of tail-anchored proteins (GET) pathway targets and inserts tail-anchored (TA) proteins into the endoplasmic reticulum (ER) membrane. A complex of two membrane proteins, Get1 and Get2, captures the TA protein from the Get3 cytoplasmic chaperone and mediates its membrane insertion.

Methods

We applied high resolution native MS based methods to recombinantly purified Get complex proteins: Get1 and Get2 membrane proteins and Get3 soluble protein. The experiments were performed on a Thermo QExactive-Plus instrument that has been modified for transmission of intact protein ions under native conditions and for trapping the ion in the source region under optional elevated activation energy (2).

Results

Overexpression and purification of Get1-Get2 membrane proteins was optimised and a full detergent screen was performed to achieve correct stoichiometry of the complex. We show that Get1-Get2 forms a dimer of dimers *in surfo*, providing the first direct evidence for their stoichiometry. With the help of conventional lipidomics approach we showed that the complex specifically interacts with a molecule of phosphatidylinositol and is destabilised by detergents that disrupt this interaction, indicating lipid-protein interactions are important within the Get1-Get2 complex. A separately purified Get3 protein was then mixed with Get1-Get2 complex and we demonstrate that Get1-Get2-Get3 interacts with a 2:2:4 stoichiometry and that each pair of

identical membrane proteins bind one Get3 dimer, providing new insights into the structure and function of the TA protein targeting complex at the ER membrane.

Conclusions

Our native MS and denaturing lipidomics data explain the role of detergents and lipids in the stability and correct stoichiometry of membrane protein translocation pathway complex. These observations for the first time explain insights into the mechanism of TA proteins translocation to ER membranes and provide a platform for further biophysical characterisation of the complex by X-ray crystallography or electron microscopy.

Novel Aspect

MS potential for understanding the fundamental biological processes, in this case, Get complex is revealed and the study opens a door for membrane protein stoichiometry related questions.

References

1. Mehmood, S. et al., *Annual Review of Physical Chemistry*, 66, 453-74 (2015).
2. Mehmood, S. et al., *Nature Chemistry*, 8, 1152-58 (2016).

WOr.01-3 CHARGE-INDEPENDENT MASS SPECTROMETRY OF SINGLE VIRUS CAPSIDS ABOVE 100MDA WITH NANOMECHANICAL RESONATORS

Christophe Masselon⁽¹⁾ - ***Sebastien Hentz***⁽²⁾ - ***Sergio Dominguez-Medina***⁽¹⁾ - ***Shawn Fostner***⁽²⁾ - ***Mohammad Halim***⁽¹⁾ - ***Martial Defoort***⁽²⁾ - ***Emeline Vernhes***⁽³⁾ - ***Pascale Boulanger***⁽³⁾

⁽¹⁾ CEA, DRF/BIG/Large Scale Biology Laboratory, Grenoble - ⁽²⁾ CEA, DRT/LETI/DCOS, Grenoble - ⁽³⁾ I2BC, Université Paris-Saclay, Gif sur Yvette

Keywords: Nanomechanical resonators, Native MS, instrumentation, polymer, virus.

Introduction

Nanomechanical resonators can measure the mass of individual particles accreting on their surface [1]. As their frequency-to-mass relationship scales inversely with the resonator's dimension, research groups have engaged in a race to measure the smallest detectable mass with increasingly small resonators [2]. However, nanomechanical resonators have best performances in the MDa to GDa mass range, where most technologies, including MS, struggle to measure masses.

Methods

We developed a novel architecture for nanoresonator-based MS of high-mass biological nanoparticles combining solution nebulization by surface acoustic wave nebulization, efficient transfer and focusing by an aerodynamic lens without relying on electromagnetic fields, and mass measurement using an array of frequency-addressed nanomechanical resonators.

Results

The system's particle transfer and focusing performances were characterized using NIST standard 45 nm polystyrene nanoparticles. Particle transfer from solution to the detector chamber reached ~7.5%, and the particle beam could be collimated to a 1.5 mm diameter with divergence better than 1°. Mass measurements of NIST nano-particles were in agreement with expectations based on polystyrene density and specified size distribution. The system was then employed for the measurement of empty and filled bacteriophage viral capsids with theoretical masses up to 105.4 MDa demonstrating its bio-analytical capabilities. Measured masses were within ~1% of their theoretical values for both empty and filled capsids. The measurement of the mass of bacteriophage T5 capsids containing the viral genome is, to the best of our knowledge, the largest defined molecular mass ever measured, and with a resolution of 75, it is also the highest reported resolution for single particle measurements above 10 MDa.

Conclusions

We have successfully developed a new nanoresonator-based MS system for high mass measurement. This system determined the mass distribution of ~30 MDa polystyrene nanoparticles with a detection efficiency 6 orders of magnitude higher than previous nanomechanical-MS systems with ion guides, and successfully performed the highest molecular mass measurement to date with less than 1 picomole of bacteriophage T5 105 MDa viral capsids.

Novel Aspect

An efficient nanoresonator-based MS system architecture allowed individual particle mass measurement without ionization in the MDa to GDa mass range.

References

1. E. Sage et al., *Nat. Commun.* 6, 6482 (2015).
2. K. Jensen, K. Kim, A. Zettl, *Nat. Nanotechnol.* 3, 533 (2008).

WOr.01-4 A CROSS-LINKING/MASS SPECTROMETRY WORKFLOW BASED ON MS-CLEAVABLE CROSS-LINKERS AND THE MEROX SOFTWARE FOR MAPPING PROTEIN-PROTEIN INTERACTIONS

Claudio Iacobucci - Michael Götze - Andrea Sinz

Martin-Luther-Universität Halle-Wittenberg, Institut für Pharmazie, Halle

Keywords: Cross-linking/Mass Spectrometry, Structural proteomics, Protein Interaction, Software, Collisional Activation

Introduction

Mass spectrometry (MS)-cleavable cross-linkers are moving into the focus of the cross-linking/MS approach for studying 3D-structure of proteins and protein assemblies. They undergo specific fragmentations under collisional activation conditions generating characteristic product ion patterns. They have an enormous potential for a fast and reliable identification of cross-linked peptides even from highly complex samples, e.g. whole cells and organisms.

Methods

We developed a robust and widely applicable workflow that allows a facile identification of cross-links for deriving spatial constraints from proteins and protein complexes. Our protocol combines the synthesis of novel cross-linkers, protein cross-linking, bottom-up proteomics analysis, and data analysis based on the in-house MeroX software [1]. Two Orbitrap mass spectrometers, Orbitrap Fusion Tribrid and Q-Exactive Plus, were employed.

Results

We designed and synthesized four novel MS-cleavable (photo)-cross-linkers. 2,2'-Azobis(2-methylpropanimidate) is an innovative cross-linking principle that after collisional activation induces free radical-initiated sequencing (FRIPS) of connected peptides in positive ion mode [2,3]. 1,1-carbonyldiimidazole is the first zero-length MS-cleavable linker for nucleophilic groups (amines and hydroxyls) in proteins [4]. 3,4-Diallylurea is a novel photo-activatable, MS-cleavable reagent to selectively target cysteines [5]. We unveiled the true reactivity of diazirine-base cross-linkers, which are in fact photo-activatable and MS-cleavable reagents to target carboxylic acids in proteins [6]. All cross-linkers were investigated to study several protein systems, such as the tumor suppressor p53, which is an intrinsically disordered protein, and the whole *E. coli* ribosome. We are currently extending our cross-linking/MS workflow to the *in-vivo* analysis of protein-protein interactions, exemplified for *Drosophila* embryos.

Conclusions

Our integrated cross-linking workflow allows to map protein 3D-structures and protein-protein interactions *in-vitro* and *in-vivo*. We have synthesized innovative, MS-cleavable cross-linking principles that target different functional groups in proteins, such as amines, hydroxyls, thiols or react non-specifically with all 20 proteinogenic amino acids. The MeroX software is able to analyze cross-linked products in an automated fashion.

Novel Aspect

We present an integrated workflow based on novel MS-cleavable cross-linking principles and a fully automated analysis of cross-linked products by the MeroX software.

References

1. Götze M., Pettelkau J., Fritzsche R., Ihling C. H., Schäfer M., Sinz A. J. Am. Soc. Mass Spectrom., 26, 83-97 (2015).
2. Iacobucci C., Hage C., Schäfer M., Sinz A. J. Am. Soc. Mass Spectrom., 28, 2039-2053 (2017).
3. Iacobucci C., Schäfer M., Sinz A., Mass Spectrom. Rev., DOI: 10.1002/mas.21568 (2018).
4. Hage C., Iacobucci C., Rehkamp A., Arlt C., Sinz A. Angew. Chem. Int. Ed., 56, 14551-14555 (2017).
5. Iacobucci C., Piotrowski C., Rehkamp A., Ihling C. H., Sinz A., J. Am. Soc. Mass Spectrom., DOI: 10.1007/s13361-018-1952-8 (2018).
6. Iacobucci C., Götze M., Piotrowski C., Arlt C., Rehkamp A., Ihling C. H., Hage C., Sinz A, Anal. Chem., 90, 2805-2809 (2018).

WOr.01-5 CAN WE USE HDX-MS TO MODEL PROTEIN CONFORMATIONS?

Antoni Borysik

King's College London, King's College London, London

Keywords: Hydrogen deuterium exchange, protein structure determination

Introduction

Protein structures are commonly utilized to interpret hydrogen deuterium exchange mass spectrometry (HDX-MS) data not vice versa. But the ability to model structures by HDX-MS would exploit the many advantages of the technique particularly with regard to speed and sensitivity. The modelling of protein folds and homo/hetero protein assemblies by HDX-MS is described and the utility and limitations of the approach for different protein systems is discussed.

Methods

HDX-MS data of proteins and their assemblies were obtained by continuous labelling HDX-MS. Computational approaches were used to simulate HDX-MS data directly from protein structures or docking outputs. A range of structural analysis/evaluation techniques were employed to quantify the efficacy of the approach to model protein folds and assemblies including hierarchical clustering and binary structural classification.

Results

The HDX-MS patterns of proteins simulated directly from their atomic structures are compared to experimental data obtained in-house. Utilizing a simple expression to estimate HDX protection factors we show that it is possible to simulate protein HDX-MS patterns to sufficient accuracy to permit the modelling of protein conformations. Protein folds and hetero-assemblies can be modelled to high accuracy but the modelling of homo-assemblies is significantly more challenging due to the inability to prepare difference plots for these systems. Surprisingly, even

poorly simulated HDX-MS data is shown to be efficient at identifying native conformations from a high background of decoy states. And in the case of homo-proteins the diagnostic ability of the HDX-MS simulations is shown to scale with the number of subunit interactions in the protein complex. Taken together these data add to our understanding of the use of HDX-MS for structural evaluation and provide an important foundation on which future developments in the area can be built.

Conclusions

HDX-MS has potential in protein modelling with applications in areas such as ab initio protein folding and protein-protein docking. [1] HDX-MS data can be simulated directly from protein structures significantly enhancing the throughput and applicability for quantitative ranking of modelling outputs. The present work sets a new precedent for the analysis of HDX-MS outputs and establishes greater synergy of the method with molecular modelling techniques.

Novel Aspect

The ability to utilize HDX-MS data to model protein structures represents a new way of working with the technique.

References

1. Borysik A. J., *Angew. Chem. Int. Ed. Engl.*, 56, 9396-9399 (2017)

WOr.02 - JMS Award Symposium

Chairs: Richard Caprioli, John Langley

WOr.02-1 ANION FORMATION IN MALDI DEPENDS ON MATRIX AND TARGET PLATE MATERIAL CHOICE

Guido Paul Zeegers - Renato Zenobi

ETH Zürich, Chemistry and Applied Biosciences, Zürich

Keywords: Target plate material, MALDI, anion formation

Introduction

In MALDI mass spectrometry, initial charge separation and charge transfer processes remain elusive. To improve analyte anion signal intensity, enhance sensitivity, and gain more insight in the roles of the target plate material and matrix in MALDI anion formation, a systematic study was carried out. The influence of 21 different metals and alloys substrates, and a range of matrices, in the absence and presence of diphenyl phosphate (DPP), on anion yield, was studied.

Methods

Molecular anion signal intensities were systematically analyzed and compared by varying laser fluence for substrates with varying melting points, thermal conductivity, electrical resistivity and work function. A comparison of DPP molecular anion signal intensity was made, based on the matrix UV-absorption (λ : 337.1 nm), electron affinity and crystallization. As many factors as possible were considered to explain observed anion signal intensity differences.

Results

Metals with oxide layers that, upon surface ablation, show an intrinsically high anion formation in absence of a deposited sample, provide the highest DPP anion signal intensity. Using Ag as target plate material generally increases the analyte signal intensity 10-fold relative to the conventionally used stainless steel. Metals with high melting points (such as Cr, W, Mo, Ta) tend to underperform. For most metals, the homogeneously deposited fullerene-C₆₀ performed best

as a matrix, though anthraquinone and 9-cyanoanthracene performed best on Ag at the highest tested laser fluence.

Conclusions

We conclude that when matrix and target plate material are chosen carefully, this can significantly enhance the ion signal intensity of an added analyte in MALDI. Hence, the commonly used stainless-steel target plates should be reconsidered. The target plate material acting as a catalyst, resulting in different anion signal intensities, should be considered, along with the potential of MALDI to predict catalysts for (photo-)chemical reactions.

Novel Aspect

Target plate material strongly affects the outcome of the ionization process and a proper choice enhances ion yields significantly.

WOr.02-2 **COMBINING BIO-CONJUGATION AND MATRIX-ASSISTED LASER DESORPTION IONISATION MASS SPECTROMETRY FOR BIOMOLECULE DETECTION AND LOCALISATION IN TISSUES**

Katherine Stevens⁽¹⁾ - **Parul Mittal**⁽²⁾ - **Arun Everest-Dass**⁽³⁾ - **Andrew Abell**⁽¹⁾ - **Tara Pukala**⁽¹⁾

⁽¹⁾ University of Adelaide, Physical Sciences, Adelaide - ⁽²⁾ University of Adelaide, Biological Sciences, Adelaide – ⁽³⁾ Griffith University, Institute for Glycomics, Southport

Keywords: MALDI, imaging, glycoprotein, cross-linking, bioconjugation

Introduction

Membrane proteins are central nodes in the networks that allow cells to receive and transmit signals between themselves and their environment. These proteins comprise 30% of the proteome and their diversity is enhanced by post-translational modifications, such as glycosylation [1]. These proteins are involved in a range of diseases, including various cancers and inflammatory diseases [2,3].

Methods

We have developed a matrix-assisted laser desorption ionisation (MALDI)-mass spectrometry (MS) approach that combines cross-linkers capable of targeting intact proteins within tissues, to enable their detection via MALDI imaging. This approach uses a novel photocleavable cross-linker that is released by excitation with a MALDI laser to produce a signature mass signal in the resulting spectra.

Results

We have demonstrated a simple and versatile synthetic route for developing MALDI imaging cross-linkers, which contain a photocleavable moiety that is released by excitation using a MALDI laser and produces a signature mass signal in the resulting spectra. The mass signal can be easily modified as the design allows for a one-step addition of peptides of various lengths at the site of photocleavage within the cross-linker. Furthermore, the specificity of this probe can be modified to enable detection of different target biomolecules. The practical application of these cross-linkers has been confirmed by MALDI and electrospray ionisation (ESI)-MS experiments. We have found using MALDI experiments that we are able to detect the desired mass signal, which corresponds to the release of the expected photocleavage product.

Conclusions

Based on our initial findings, we aim to apply this method for the detection of multiple glycoprotein targets simultaneously in tissue sections, as a potential new approach to cancer diagnosis. Additional ongoing research is focussed on optimising this technique for glycoprotein

detection and developing similar approaches for targeting DNA oligos and glycan-binding proteins within tissues.

Novel Aspect

This approach provides a method for detecting intact biomolecules using MALDI-MS. It has the potential for multiplexing by modifying the length of the peptide mass signal.

References

1. Moraes, Isabel, et al. *Biochimica et Biophysica Acta (BBA)-Biomembranes* 1838.1, 78-87 (2014).
2. Connelly, Margery A., et al. *Clinica Chimica Acta* 459, 177-186 (2016).
3. Banerjee, Nilotpai, and Sumi Mukhopadhyay. *Virus Disease* 27.1, 1-11 (2016).

WOr.02-3 HISTOLOGY-DIRECTED PROTEOMICS: COMBINING AUTOFLUORESCENCE GUIDED MICRO-DIGESTIONS WITH LIQUID SURFACE EXTRACTIONS

Daniel Ryan⁽¹⁾ - **Nathan Patterson**⁽²⁾ - **Jessica Moore**⁽²⁾ - **Jeffrey Spraggins**⁽³⁾ - **Richard Caprioli**⁽³⁾

⁽¹⁾ Vanderbilt University, Department of Chemistry, Nashville - ⁽²⁾ Vanderbilt University, Mass Spectrometry Resource Center, Nashville - ⁽³⁾ Vanderbilt University, Department of Biochemistry, Nashville

Keywords: Autofluorescence, LESA, Bottom-Up Proteomics

Introduction

Liquid extraction surface analysis (LESA) has proven to be an effective approach for generating spatially-targeted, proteomics data from tissue.^{1,2} Current approaches enzymatically digest the whole tissue surface prior to extraction. When foci are smaller than the LESA spot size, signal from the area of interest is diluted. Herein, we demonstrate proteomic analyses of foci using histology-directed micro-digestions guided by autofluorescence (AF) microscopy.

Methods

Rat kidney tissue was sectioned and mounted onto microscope slides. AF (fluorescence of endogenous molecules) microscopy images of a section was used to guide trypsin spotting onto foci of interest using a piezoelectric spotter. LESA samples from the individually digested foci were pooled and analyzed by LC-MS/MS on a Q Exactive Orbitrap mass spectrometer. Adjacent sections were digested using robotically sprayed trypsin to compare to traditional LESA approaches.

Results

Proof-of-concept experiments were completed using LESA to extract from undigested control sections, sections that were robotically sprayed with a homogenous coating of trypsin, and sections with trypsin spotted on foci guided by autofluorescence microscopy. The control experiments resulted in ~10 protein IDs per extract. Robotically sprayed tissue sections contained on average 250 protein identifications per mm² sampled by LESA. The robotic spraying approach also suffers from some analyte delocalization during the spray process, as monitored by AF microscopy of the tissue following digestion. This approach is poorly suited for targeting foci smaller than the size of the LESA extraction. LESA samples from micro-digested kidney using AF-guided trypsin spotting provided upwards of 1500 protein ID's per mm² sampled by LESA, a 5-fold increase compared to the robotic spray method. The diameter of the trypsin spots approach 100 μm at its smallest, providing a 3-4 fold improvement of spatial resolution when compared to extracting via LESA alone.

Conclusions

LESA is limited to relatively large droplet diameters and is inadequate for studying small, discrete foci (< 400 µm) often of interest in biological tissues. The combination of tryptic micro-digestions guided by AF microscopy and LESAs circumvents LESAs spot size limitations by only producing digested peptides from the targeted foci.

Novel Aspect

Autofluorescence-guided tryptic micro-digestions combined with LESAs allows for targeted proteomic analysis of small biological foci in tissue.

References

1. Lamont, L.; Baumert, M.; Potocnik, N. O.; Allen, M.; Vreeken, R.; Heeren, R. M. A.; Porta, T., *Analytical Chemistry* **2017**, *89* (20), 11143-11150.
2. Sarsby, J.; Martin, N. J.; Lalor, P. F.; Bunch, J.; Cooper, H. J., *Journal of the American Society for Mass Spectrometry* **2014**, *25* (11), 1953-1961.

WOr.02-4 DUAL MICROPROBE CE-ESI-HRMS FOR THE CHARACTERIZATION OF PROTEINS AND METABOLITES IN SINGLE EMBRYONIC CELLS OF LIVE FROG EMBRYOS

Camille Lombard-Banek⁽¹⁾ - Rosemary Onjiko⁽¹⁾ - Sally Moody⁽²⁾ - Peter Nemes⁽¹⁾

⁽¹⁾ University of Maryland, Chemistry & Biochemistry, College Park - ⁽²⁾ The George Washington University, Anatomy and Regenerative Biology, Washington

Keywords: Single-cell, CE-ESI-HRMS, Systems biology, Proteomics, Metabolomics

Introduction

Establishment of cell-to-cell heterogeneity is critical to the normal formation of different tissues in vertebrates. A system biology understanding of this heterogeneity with single cell resolution could provide important insights on how cell fate is established. We thereby developed a dual-microprobe sampling approach hyphenated to CE-ESI-HRMS to measure proteins and metabolites from the same identified single embryonic cells in the 8-cell frog embryo.

Methods

A 10 nL volume was withdrawn from identified cells of the 8-cell *X. laevis* embryo using capillary microsampling. Metabolites were extracted in 40% acetonitrile/40% methanol solution. Proteins were delivered into 50 mM ammonium bicarbonate and trypsin-digested. All extracts were measured using our custom-built CE-ESI-HRMS platforms. Metabolites were quantified by using peak-area under the curve as a proxy of relative concentration. Proteins were quantified using MaxLFQ approach integrated in MaxQuant.

Results

To advance the systems biology understanding of cell heterogeneity during embryonic development, we developed dual-microprobe sampling to collect two different portions of the cytoplasm from the same cell to measure metabolites and proteins. We used this approach on the D1 and V1 cells in the same 8-cell embryo. These cells are situated across the dorsal/ventral axis of the embryo. The approach is based on our previous microprobe-CE-ESI-HRMS protocol, whereby an ~10 nL volume of the cell cytoplasm is extracted, leaving the embryo essentially intact.¹⁻³ We consistently quantified 335 proteins and 150 metabolite features. The different cells separated using PLS-DA analysis based on protein and metabolites profiles. Overall, 6 proteins and 5 metabolite features were differentially expressed. Using these quantitative data, we performed pathway enrichment analysis. We identified the arginine and proline metabolism

pathway as one of the most statistically significantly enriched, which suggests its importance in early embryonic development.

Conclusions

We have developed a dual sampling approach enabling the analysis of proteins and metabolites from the same cell. The collected data represent the first example of proteomic and metabolomic analyses from the same single-cell in a live embryo, providing new information on cell heterogeneity in the 8-cell frog embryo. Moreover, the technology promises to deepen our understanding of the mechanisms involved in embryo patterning during embryonic development.

Novel Aspect

Dual microprobe-CE-ESI-HRMS enabled proteomic and metabolomic analyses in the same single-cell of live frog embryos.

References

1. Onjiko, R. M.; Portero, E. P.; Moody, S. A.; Nemes, P., *Anal. Chem.*, 89, 7069-7076 (2017)
2. Onjiko, R. M.; Plotnick, D. O.; Moody, S. A.; Nemes, P., *Anal. Methods*, 9, 4964-4970 (2017)
3. Lombard-Banek, C.; Moody, S. A.; Manzini, M. C.; Nemes, P., *Submitted*

WOr.02-5 WASTEWATER-BASED EPIDEMIOLOGY: USE OF ADULTERANTS TO ESTIMATE COCAINE CONSUMPTION PATTERNS IN BRASÍLIA, BRAZIL BY LC-QTOF-MS

Katyeny Manuela Da Silva⁽¹⁾ - **Adriano Otávio Maldaner**⁽²⁾ - **Fernando Fabríz Sodré**⁽¹⁾

⁽¹⁾ University of Brasília, Chemistry Institute/Faculty of Technology, Brasília - ⁽²⁾ Federal Police, Federal Police, Brasília

Keywords: waste water epidemiology, drugs of abuse, forensic toxicology, cocaine, adulterants

Introduction

Illicit drugs are a major public safety problem in Brazil. Findings of the Brazilian Federal Police concerning drug seizures show that levamisole (LEV) is the main adulterant in cocaine hydrochloride (mean of 13.2%), while phenacetin (PHE) is the mainly found in freebase cocaine – crack (mean of 15.2%).^{1,2} Thus, strategies based on quantification of LEV and PHE, together with the major cocaine metabolite, benzoylecgonine (BE), in wastewater may contribute to identify different forms of the consumed cocaine.³

Methods

Composite (24 h) raw sewage samples were collected at the entry of the two main wastewater treatment plants (NW-WWTP and SW-WWTP) in the Brazilian Capital, Brasília. Aliquots of 50 mL (pH 2) were passed through solid phase extraction cartridges (StrataX) previously conditioned with MeOH:ACN (60:40) and H₂O (pH 2).³ Analytes were recovered with 6 mL of MeOH, the extracts were spiked with isotopically labeled standards (BE-d₃ and COC-d₃) and concentrated to 0.2 mL. Analyses were carried out using a liquid chromatography system coupled to a quadrupole-time-of-flight mass spectrometry (LC-QTOF).

Results

Results revealed a cocaine consumption varying from 750±120 to 2470±410 mg/day/1000 inhab in weekdays and 2340±1190 to 4580±750 mg/day/1000 inhab in weekends. In terms of load, considering cocaine purity in Brazilian seizures (49 ± 29 %)², we estimate an annual consumption of 2 tons of cocaine, in the region. PHE/LEV ratios found in seized cocaine street samples varies from 0.18 in hydrochloride cocaine to 15.2 in freebase samples, indicating the prevalence of crack consumption as PHE/LEV ratios becomes higher. Relationships between both adulterants

in the samples evidenced a higher use of crack cocaine in the NW-WWTP region (PHE/LEV = 3.0 ± 1.5) compared to SW-WWTP. The PHE/LEV ratio was the same in the SW-WWTP region (PHE/LEV = 1.5 ± 1.0) when compared to 2014 data,³ while the ratio for NW-WWTP region doubled suggesting an increase of crack users in this region.

Conclusions

The LC-QTOF-MS system allows the identification based on accurate mass measurements and quantification with multiple reaction monitoring. Regarding the forensic field, it is very important to know how a drug is consumed. The use of adulterants, such as PHE and LEV may be useful to provide novel information concerning cocaine consumption patterns and polyuse of drugs.

Novel Aspect

Alternative approach to differentiate crack and hydrochloride cocaine consumption based on the quantification of adulterants in wastewater samples.

References

1. Botelho E. D., Cunha R. B., Campos A. F. C., Maldaner A. O., *J. Braz. Chem. Soc.* **2014**, *25*, 611.
2. Maldaner A. O., Botelho E. D., Zacca J. J., Melo R. C. A., Costa J. L., Zancanaro I., Oliveira C. S. L., Kasakoff L. B., Paixão T. R. L. C., *J. Braz. Chem. Soc.* **2016**, *27*, 719.
3. Sodr  F. F., Souza G. B., Feitosa R. S., Pereira C. E. B., Maldaner A. O., *J. Braz. Chem. Soc.* **2017**, *28*, 2146.

WOr.03 - Reaction mechanisms and elusive intermediates

Chair: Scott A. McLuckey

WOr.03-1 **Keynote: ELUCIDATING THE CHEMISTRY OF ELUSIVE GAS PHASE RADICALS USING DISTONIC RADICAL IONS: TRYING TO GET A GOOD DEAL FROM A FAUSTIAN BARGAIN**

Stephen Blanksby - David L. Marshall - Berwyck L. J. Poad - Adam Trevitt

Central Analytical Research Facility (CARF), Institute for Future Environments, Queensland University of Technology

Keywords: ion-trap mass spectrometry; distonic radical ions; ion spectroscopy; thermochemistry

Introduction

Distonic radical ions can be used to explore the chemistry of transient radicals in the gas phase. While this approach facilitates the experimental observation of elusive radicals, it also perturbs their native reactivity. Direct thermodynamic and spectroscopic measurements of the influence of charge on radical stability within distonic radical ions have been undertaken and highlight perturbations over significant intramolecular separations.

Methods

Radical ions have been synthesized in the gas phase using regioselective approaches based on collision-induced dissociation and photodissociation on modified ion trap mass spectrometers. The kinetics of ion-molecule reactions have been measured by seeding low concentrations of reagent gas in the helium buffer gas supplied to the ion trap.

Results

Pioneering research by Kentt maa has demonstrated the power of the distonic ion approach but has also highlighted the Faustian bargain that is entered into with the introduction of the

charged moiety to the radical [1]. That is, while the presence of the charge facilitates the experimental observation and characterization of elusive radicals it can also perturb their native energetics. In our laboratory a homologous series of distonic ions incorporating nitroxyl and peroxy radicals have been regioselectively prepared as both positively and negatively charged ions. The influence of the charge polarity and the distance between the charge and radical sites were systematically explored using a combination of ion spectroscopy; thermochemical measurements (using the kinetic method); and ion-molecule reactions. Key findings from these investigations reveal that radical moieties are stabilised in distonic anions and destabilised in analogous cations, with respect to their corresponding neutrals. Significantly, these effects could be measured over intramolecular distances of up to 10 Å.

Conclusions

The thermodynamic impact of charged groups of radical stability provides the critical insight required to extrapolate distonic ion behaviours to their neutral archetypes. Moreover, the long range nature of charge-radical effects suggests that these phenomena are likely playing a significant role in the dissociation chemistries of the radical ion intermediates formed in contemporary ion dissociation technologies (e.g., ETD, ECD and RDD).

Novel Aspect

Direct measurement of charge-radical interaction energetics in distonic ions

References

1. Petzhold C.J., Nelson, E.D., Lardin H.A., Kenttamaa H., J. Phys. Chem. A, 106, 9767 (2002)

WOr.03-2 ACETONE AS PHOTO-INITIATOR FOR DISULFIDE REDUCTION – APPLICATION TO STRUCTURAL ANALYSIS OF DISULFIDE RICH PEPTIDES AND PROTEINS

Yu Xia⁽¹⁾ - Sarju Adhikari⁽²⁾ - Xiaoyue Yang⁽¹⁾

⁽¹⁾ Tsinghua University/Chemistry, Beijing - ⁽²⁾ Purdue University/Chemistry, West Lafayette

Keywords: Photo-Initiator, Radical Reaction, Disulfide Peptide, Tandem Mass Spectrometry, Proteomics

Introduction

Pinpointing disulfide linkage pattern is necessary for characterizing proteins containing disulfide bonds. Often carefully controlled enzymatic digestion and multi-step partial disulfide reduction are needed before MS analysis to achieve confident assignment. Herein, we describe a new method based on coupling radical initiated disulfide bond cleavage with tandem mass spectrometry for rapid and detailed structural analysis of disulfide-rich peptides.

Methods

A flow microreactor consisted of a 254 nm low-pressure mercury lamp (BHK, Inc., Ontario, CA) in parallel to UV-transparent fused silica capillary (100 µm i.d., 375 µm o.d) at 0.5 cm distance. The peptide/protein solution was infused to the reactor by a syringe pump and subjected to photochemical reactions before online electrospray ionization (ESI)-MS analysis. All data were collected on QTRAP 4500 and X500R QTOF mass spectrometers (SCIEX, Concord, ON, Canada).

Results

Acetone (1% volume ratio) was used as photo-initiator (PI) to induce subsequent radical reactions. Upon 254 nm UV irradiation, acetone underwent Norrish type I cleavage, producing acetyl radical and methyl radical[1]. The presence of alkyl alcohol as co-solvent in the peptide solution converted the primary radicals into hydroxyalkyl radicals, which subsequently cleaved

disulfide bond via substitution [2]. Disulfide reduction was the dominant reaction channel when using isopropyl alcohol, which completed in less than 5 s with >95% yield. More importantly, partial reduction of disulfide bonds could be achieved simply by shortening the reaction time. For structural characterization of disulfide-rich peptides/proteins, sequence information and disulfide bond connecting pattern could be obtained from coupling complete or partial disulfide reduction with MS/MS. The analytical utility of this method was demonstrated with peptides/proteins containing 2 to 4 disulfide bonds.

Conclusions

Disulfides are highly reactive toward radical attack; however, uncontrolled radical reactions often lead to disulfide bond scrambling, which is undesirable for disulfide linkage pattern determination. In this work we demonstrated that controlled disulfide cleavage via hydroxyalkyl radical could be achieved. Furthermore, this reaction can be directly coupled with ESI-MS/MS analysis, allowing characterization of disulfide rich peptides.

Novel Aspect

Photochemical reduction of disulfide bond in peptides and proteins and its online coupling with ESI-MS/MS for disulfide linkage pattern determination.

References:

1. Norrish, R. G. W., Bamford, C. H., *Nature*, 140, 195–6 (1937)
2. Stinson C.A., Xia Y., *J. Am. Soc. Mass Spectrom.*, 25, 1192-1201 (2014)

WOr.03-3 KINETICS OF A CLICK REACTION IN CHARGED MICRODROPLETS WITH TIME RESOLVED REACTIVE LASER ABLATION ELECTROSPRAY IONIZATION MASS SPECTROMETRY

Fred van Geenen - Satesh Gangarapu - Maurice Franssen - Michel W.F. Nielen- Han Zuilhof

Wageningen University, Laboratory of Organic Chemistry, Wageningen

Keywords: cycloaddition, kinetics, laser ablation electrospray ionization, time resolved mass spectrometry

Introduction

Time resolved mass spectrometry (TRMS) provides insights into fundamental processes in (bio)molecular reactions, such as the direct detection of short-lived intermediates and the acquisition of kinetic data. Laser ablation based injection coupled to electrospray ionization (LAESI) can be used as a novel TRMS method to perform reactions in charged microdroplets, *e.g.* to accelerate synthesis or to study the chemistry of life.

Methods

A Protea Biosciences LAESI DP-1000 system was electronically connected to a Waters Synapt G2-S time of flight mass spectrometer (MS). The distance between LAESI ion source and MS inlet was variable between 1 and 100 cm. Transport of ions from the ionization source to the MS inlet was assisted by a nebulizer gas via a grounded metal tube. Temporal resolution was obtained by the difference in time between the pulsed laser pulse (injection) and MS detection of product.

Results

We demonstrate a reactive LAESI-TRMS method with the study of a cycloaddition reaction between dipyrindyl-tetrazine and an amine derivative of bicyclononyne (BCN-amine) in charged microdroplets. For this reaction to occur a substantial increase of reaction time is required. This increase was achieved using charged microdroplet chemistry in combination with the insertion

of a reaction tube at ambient conditions between the LAESI source and the MS inlet, resulting in near complete conversions. Kinetic data was obtained with the use of temporal resolution, which was demonstrated in the 1 to 500 ms time range. Interestingly, no drop was detected in product absolute intensity over a distance up to 100 cm (0.5 s). Furthermore, this approach showed excellent instrumental stability, improved sensitivity and was demonstrated with other click reactions. Finally, reactive LAESI mass spectrometry imaging (MSI) experiments were conducted on an image containing dipyrindyl-tetrazine solution.

Conclusions

We have developed a novel reactive LAESI method for the determination of reaction kinetics in charged microdroplets and showed improved sensitivity following online reactions in LAESI analysis. Additionally, we demonstrated our methodology in a reactive LAESI-MSI experiment.

Novel Aspect

Coupling of laser ablation injection with ESI-MS to study reactions in microdroplets with temporal resolution, and demonstrating reactive LAESI-MS(I) for the first time.

WOr.03-4 AMINO-ACIDS OXIDATION: A COMBINED STUDY OF CYSTEINE OXO-FORMS BY MASS SPECTROMETRY, IRMPD SPECTROSCOPY AND SIMULATIONS

Debora Scuderi⁽¹⁾ - ***Enrico Bodo***⁽²⁾ - ***Simonetta Fornarini***⁽³⁾ - ***Barbara Chiavarino***⁽³⁾ - ***Maria Elisa Crestoni***⁽³⁾ - ***Veronica Macaluo***⁽⁴⁾ - ***Riccardo Spezia***⁽⁵⁾

⁽¹⁾ Laboratoire de Chimie Physique, UMR 8000, Université Paris Sud, Orsay - ⁽²⁾ Dipartimento di Chimica, Università degli Studi di Roma La Sapienza, Rome - ⁽³⁾ Dipartimento di Chimica e Tecnologie del Farmaco, Università degli Studi di Roma La Sapienza, Rome - LAMBE, ⁽⁴⁾ Université d'Evry, CNRS, CEA, Université Paris-Saclay, Evry - ⁽⁵⁾ Laboratoire de Chimie Théorique, UMR 7616 CNRS, Sorbonne Université, Paris

Keywords: IRMPD spectroscopy; cysteine oxo-forms; collisional induced dissociation; ab initio molecular dynamics simulations; chemical dynamics simulations

Introduction

Cysteine residues in sulfenic (-SOH), sulfinic (-SO₂H) and sulfonic (-SO₃H) forms are key intermediates in the redox-switch chemistry of proteins [1]. Although sulfenic acids have been identified in several biochemical frameworks, they are transient and not isolable in small molecular systems. Low-energy collision-induced dissociation of L-Cysteine sulfate affords L-cysteine sulfenic acid.

Methods

Chemical dynamics simulations [2] have been employed to shed light on fragmentation mechanisms of L-Cysteine sulfate. A combined approach based on infrared multiple photon dissociation spectroscopy, calculations of IR frequencies and finite temperature ab initio molecular dynamics simulations has been employed to characterize the gas-phase structures of deprotonated cysteine oxo-forms.

Results

IRMPD spectra of cysteine oxo-forms have been recorded in the highly informative mid-IR spectral range employing a Paul ion-trap mass spectrometer coupled to the IR beamline of the tunable free-electron laser at the CLIO facility in Orsay [3]. These ions show different structural motifs owing to the preferential binding of the proton to either the carboxylate or the sulfur-containing group. Due to the decreasing basicity of the sulfenic, sulfinic and sulfonic terminals, the proton is bound to SO⁻ in [cysSO]⁻, migrates on the carboxylate in [cysSO₃]⁻, while it turns out to be shared in [cysSO₂]⁻. Evidence is gathered that a mixture of close-lying low-energy

conformers is sampled for each cysteine oxo-form in a Paul ion-trap at room temperature [4]. An exotic roaming mechanism has been evidenced in the formation of L-cysteine sulfenic acid with other competing pathways.

Conclusions

The spectroscopic signatures of the oxidative modifications enable the assignment of the deprotonation site and conformational landscape of the sampled [CysSO_x] (x= 1,2,3) oxo-forms and may provide diagnostic probes to reveal critical intermediates in protein redox regulation.

Novel Aspect

The type of mass spectrometry instrument employed for collisional activation of L-Cysteine sulfate resulted to be “pathway” selective.

References

1. Reddie, K.G; Carroll K.S., *Curr Opin Chem Biol*, 12, 746-754 (2008)
2. Spezia, R.; Salpin, J.-Y.; Gaigeot, M.-P.; Hase, W. L.; Song, K. *J.Phys.Chem.A*, 113, 13853-13862 (2009)
3. Lemaire, J.; Boissel, P.; Heninger, M.; Mauclaire, G.; Bellec, G.; Mestdagh, H.; Simon, A.; Le Caer, S.; Ortega, J.M. ; Glotin, F. ; Maitre, P. ; *Phys. Rev. Lett.* 89, 3-6 (2002)
4. Scuderi, D.; Bodo, E.; Chiavarino, B.; Fornarini, S.; Crestoni, M. E. *Chem. Eur. J.*, 22, 17239–17250 (2016)

WOr.03-5 LOSS OF WATER FROM POLYGLYCINES AND THEIR [b_n]⁺ IONS; IMIDAZOLONES AND OXAZOLONES

K. H. Brian Lam⁽¹⁾ - **Justin Kai-Chi Lau**⁽¹⁾ - **Cheuk-Kuen Lai**⁽¹⁾ - **Ivan K. Chu**⁽²⁾ - **Jos Oomens**⁽³⁾ - **Alan C. Hopkinson**⁽¹⁾ - **K. W. Michael Siu**⁽⁴⁾

⁽¹⁾ York University, Department of Chemistry and Centre for Research in Mass Spectrometry, Toronto - ⁽²⁾ University of Hong Kong, Department of Chemistry, Hong Kong - ⁽³⁾ Radboud University, Institute for Molecules and Materials, FELIX Laboratory, Toernooiveld - ⁽⁴⁾ University of Windsor, Department of Chemistry and Biochemistry, Windsor

Keywords: Loss of water; imidazolone; oxazolone; [b_n]⁺ ions; interconversion between isomers

Introduction

Loss of water in the fragmentation of protonated tetraglycine is mainly from the peptide backbone, resulting in [b₄]⁺ ions with imidazolone structures. These [b₄]⁺ ions have identical CID spectra indicating isomerization prior to dissociation. Here, the fragmentations of larger protonated polyglycines [Gly_n + H]⁺ and their [b_n]⁺ ions (n=5–10) are examined. Isomerization of the [b_n]⁺ ions and their fragmentation mechanisms are studied.

Methods

All experiments were carried out using an Orbitrap Elite Hybrid Ion Trap-Orbitrap mass spectrometer. Isotopically labeled peptides were synthesized and used to determine the structures of the fragment ions. DFT calculations at B3LYP/6-311++G(d,p) level provided reaction energetics and insight into the dissociation mechanisms.

Results

¹⁸O-labeling experiment shows that loss of water in the fragmentation of protonated pentaglycines are predominantly from the first and second amide groups. Unlike the [b₄]⁺ ion derived from tetraglycines, the [b₅]⁺ ions show similar CID spectra but with different relative abundances. This suggests that the [b₅]⁺ ions undergo fragmentations more easily than isomerize. This conclusion is supported by DFT calculations that give the energy barrier for the

isomerization to be higher than that to fragmentation. Hexaglycine loses water from the first three amide groups. The $[b_6]^+$ ions derived from the water loss from the first and second amide oxygens shows similar fragmentation patterns to the corresponding $[b_5]^+$ ions. However, further water loss from the peptide backbone is found from the $[b_6]^+$ ion that was formed when water was lost initially from the third amide group.

Conclusions

The $[b_4]^+$ ions from tetraglycine possess imidazolone structures and undergo isomerization prior to fragmentation. However, the larger $[b_n]^+$ ions have different fragmentation patterns indicating that the barriers to isomerization is high. These $[b_n]^+$ ions still feature imidazolone structures, and there is no evidence for oxazolone structures.

Novel Aspect

Loss of water from protonated polyglycines are mainly from the backbone. With the increase of the peptide chain length, multiple water loss is observed.

WOr.04 – Pharmaceuticals

Chairs: Christine Eckers, Gianluca Bartolucci

WOr.04-1 **Keynote: PHARMACEUTICAL INDUSTRY CHALLENGES WITH MASS SPECTROMETRY SOLUTIONS: FROM MICROSAMPLING TO MONOCLONAL ANTIBODIES**

Walter Korfmacher

Keywords: Pharmaceuticals, Biologicals, Microsampling, HRMS, LC-MS

Introduction

The pharmaceutical industry is undergoing many changes. One challenge is to follow the 3R's by reducing the number of laboratory animals needed in the new drug discovery process. One solution is to take smaller blood samples (microsampling) so fewer laboratory animals are needed [1]. Another challenge is the shift toward biologicals. Biologicals are more difficult to assay than small molecules.

Methods

For blood microsampling, either capillary microsampling (CMS) tubes or a volumetric absorptive microsampling (VAMS) device was used. For mass spectrometry imaging (MSI), various high resolution mass spectrometer (HRMS) systems are discussed. For biological assays, various MS systems are described including triple quadrupole systems and HRMS systems.

Results

In one study, we used the 8- μ L CMS blood sampling method [2]. In another study, we demonstrated the utility of both blood CMS and VAMS for getting samples from rat PK studies [3]. MSI has now become an important tool for understanding the distribution of a drug as well as its metabolites in laboratory animals. MSI can also be utilized for identifying disease and toxicity biomarkers as part of new drug discovery. Various strategies have been developed to measure therapeutic biological compounds, including monoclonal antibodies, in serum from mouse or monkey PK studies. One strategy is the "surrogate peptide approach", where the large biological compound is cut into smaller pieces (typically by using trypsin digestion). Then one or more of the peptides is selected to serve as the analyte that is quantified (typically by LC-

MS/MS) in order to then calculate the concentration of the large biological molecule. Another approach is to use LC-HRMS to measure the intact biological compound directly.

Conclusions

Overall, these studies provide additional data that blood CMS is a valuable approach for serial blood sampling from discovery rat PK studies. Mass spectrometry continues to be a key system that is used throughout all stages of the pharmaceutical industry research paradigm, from early drug discovery to clinical studies. As MS systems continue to improve, then new applications become possible.

Novel Aspect

The talk will describe how MS has been utilized to provide solutions for the pharmaceutical industry during new drug discovery and development.

References

1. Chapman K., Chivers S., Gliddon D. et al., *Drug Discovery Today* 19 (5), 528-532 (2014).
2. Korfmacher W., Fitzgerald M., Luo Y. et al., *Bioanalysis* 7(4), 449-461 (2015).
3. Luo Y., Korfmacher W., Ho S. et al., *Bioanalysis*, 7(18), 2345-59 (2015).

WOr.04-2 **IMPACT OF CENTROIDING AND DIFFERENCES IN DATA PROCESSING ON QUANTIFICATION USING HIGH-RESOLUTION MASS SPECTROMETRY**

Liesbeth Vereyken - Lieve Dillen - Rob J. Vreeken - Filip Cuyckens

Janssen R&D, Discovery Sciences, Beerse

Keywords: High-Resolution MS, quantification, data processing, mass extraction window

Introduction

The improved performance of High-Resolution MS (HRMS) systems triggers more (bio)analysts to switch from triple quad MS to HRMS for quantitative analysis. The impact of differences in data handling by different software packages, the impact of centroid or continuum data as well as the use of theoretical or measured accurate mass on the selection of the optimal narrow mass extraction window was studied.

Methods

Calibration curves of 8 pharmaceutical drugs in plasma were analyzed in 5-fold on a Synapt G2-S QToF system. Data were acquired in continuum mode and processed in four different ways: as centroid data, as individual profile data points and as full profile data using 2D and 3D processing. One compound was analyzed on 2 QTofs and one Orbitrap system and processed in six different software packages.

Results

Centroid data provided advanced selectivity resulting in cleaner XICs for qualitative analyses but with a limited to no effect on quantitative performance (LOQ, precision, linearity). The optimal mass extraction window (MEW) selection showed to be solely defined by the mass accuracy of the system used. For profile data, the improved selectivity provided by narrower MEWs is counterbalanced by a loss in MS peak area affecting signal-to-noise. Therefore, the full width half maximum (FWHM) of the peak showed to be the ideal starting point for MEW selection in addition to mass accuracy that becomes gradually more important with increasing MS resolving power. 3D peak processing can overcome the challenges of MEW selection but comes with different challenges such as the difficulty in reviewing and adjusting bad peak integrations.

Conclusions

The impact of centroid or profile data and different software approaches on the extraction of MS data to generate XICs were investigated. While there is an effect on the optimal MEW selection for which guidelines are given, the effect on quantitative performance was minimal.

Novel Aspect

The impact of differences in data treatment by software packages. Optimization of processing parameters for quantitative HRMS performance

WOr.04-3 LC-ICP-MS AS AN ALTERNATIVE TO RADIOACTIVITY IN PHARMACOLOGY

Christine Enjalbal⁽¹⁾ - ***Emmanuelle Cordeau***⁽¹⁾ - ***Gilles Subra***⁽¹⁾ - ***Didier Gagne***⁽¹⁾ - ***Nolween Prache***⁽¹⁾ - ***Carine Arnaudguilhem***⁽²⁾ - ***Brice Bouyssiere***⁽²⁾

⁽¹⁾ IBMM, University of Montpellier, Montpellier - ⁽²⁾ IPREM, UPPA, Pau

Keywords: ICP-MS; Peptide; Quantification; Pharmacology; Labeling

Introduction

Measuring the affinity of a ligand for its target is an issue of central importance in drug development. At present, radiolabeling is the reference method in pharmacology, allowing a highly sensitive detection of the molecule of interest but with many drawbacks (user safety, environmental issue). In this context, we developed an alternative methodology using mass spectrometry due to its intrinsic capabilities in terms of specificity and sensitivity of analysis.

Methods

The cornerstone of this analytical methodology relies on the production of very abundant signals related exclusively to the peptide ligand of interest to be quantified by a generic mass spectrometry method, applicable to all receptor/ligand systems, while demonstrating performance equal to radioactivity. Among all techniques investigated [1-3], LC-ICP-MS associated with peptide selenium labeling provided extreme efficacy and robustness of quantification.

Results

This strategy was applied to quantify peptide ligands at very low concentrations from cell cultures in pharmacological experiments (saturation and competitive binding assays). Noteworthy, we developed a straightforward protocol that can be performed with a standard macrobore LC-ICP-MS system without any sample treatment or special instrumental set-up, while allowing efficient suppression of all mass interferences to reach the targeted sensitivity. Significantly, a quantification limit of 50 ng Se L-1 (50 ppt/63 fmol of injected Se) was achieved, the samples issued from the pharmacological assays being directly introduced into the LC-ICP-MS system. The proof of concept was illustrated on the vasopressin/AVP pharmacological model. The methodology was then successfully applied to a second ligand/receptor system (CCK-4/CCK) proving the generic aspect of the protocol. A very good correlation of all experimental data for the determination of K_i , K_d and IC_{50} by both ICP-MS and radioactivity demonstrated the reliability of the strategy.

Conclusions

We aimed to obtain a generic mass spectrometry method, applicable to all receptor/ligand systems, while demonstrating performance equal to radioactivity in order to convince pharmacologists to modify their practices by the eradication of the use of radioactive substances. The gathered results with two different ligand/receptor systems indicate clearly that the LC-ICP-MS-based developed method represents a viable alternative to radiolabeling.

Novel Aspect

LC-ICP-MS equaled radioactivity for peptide quantification from cell cultures which should convince pharmacologists to modify their practices.

References

1. Paramelle D. et al., *Angew. Chem. Int. Ed.* 49, 8240 (2010)
2. Rossato M. et al., *Chem. Med. Chem.* 11, 2582 (2016)
3. Cordeau E. et al., *Plos One*, 11(6): e0157943 (2016).

WOr.04-4 SUB-PPM DETECTION AND IDENTIFICATION OF IN-VIVO METABOLITES OF AN ANTIBODY-DRUG CONJUGATE BY DATA-INDEPENDENT HPLC - HIGH-RESOLUTION TANDEM MS

Antonio Triolo⁽¹⁾ - **Elisa Ciuti**⁽¹⁾ - **Chiara Carboni**⁽²⁾ - **Vera Chirulli**⁽²⁾

⁽¹⁾ Menarini Ricerche, Analytical Chemistry for Development, Firenze - ⁽²⁾ Menarini Ricerche, Pharmacokinetics and Metabolism, Pomezia

Keywords: antibody-drug conjugates, anti-tumor, metabolites, data-independent, tandem MS

Introduction

The activity of an Antibody-drug conjugate (ADC) is exerted after drug release from the ADC inside the cancer cells. Drug and its cellular metabolites are further transformed in vivo, generating a series of highly cytotoxic drug-related molecules which must be identified as a first step for their pharmacological and toxicological evaluation. Here we describe identification of the payload metabolites in serum of cynomolgus administered intravenously with an ADC.

Methods

The deproteinized sera were analyzed by LC-MS/MS with a 100 x 1.0 mm RP column, using a Q Exactive mass spectrometer operated in positive ESI. Data acquisition consisted of a full-scan MS at FWHM resolution = 70000, alternated with a series of 30 MS/MS scans, at resolution = 17500 and HCD collision energy 20 V, acquired with precursors ranging sequentially from m/z 610 to 1161 in steps of 19 m/z units, and 20 m/z units quadrupole isolation windows. Scan time was 2.8 secs.

Results

The above described method afforded full-scan sensitivity to sub-ng/mL levels, and sufficient selectivity to distinguish the metabolites signals in the presence of coeluting matrix peaks up to 3 degrees of magnitude higher in abundance. The payload-related metabolites were detected by extracting the accurate masses of selected fragment ions from the entire MS/MS data set, whereas their precursors were visually searched at the corresponding retention times in the MS data set, inside the 20 m/z-wide interval from which the fragment ions were generated.

The method selectivity was verified by analyzing blank serum, and the capability of detecting unknown payload-related compounds by spiking blank serum with 5% v/v drug incubates in human liver microsomes.

Cynomolgus sera treated with 5 and 1 mg/kg ADC showed the presence of the free drug, its cysteine conjugate, and the drug plus linker adduct. Each metabolite was characterized by the accurate masses of the precursor and its main diagnostic fragment ions.

Conclusions

High-resolution MS and data-independent analysis with sequential windows MS/MS acquisition of a series of precursors spanning an about m/z 600 mass interval, allowed successful detection and identification of the payload-related metabolites in cynomolgus serum after treatment with an ADC at 5 and 1 mg/kg dose. The method features high specificity and full-scan sensitivity in the sub-ng/mL range, with mass accuracy typically below 2 ppm.

Novel Aspect

First data-independent scanning application to the detection and identification of in vivo ADC metabolites.

WOr.04-5 NEW QC IDENTITY TESTING FOR ACELLULAR PERTUSSIS COMBINATION VACCINES USING LC-MS

Hong Zhang - Rob Maharajh - Seeven Vydelingum - Lois Yin

Sanofi Pasteur, Analytical Process and Technology, Toronto

Keywords: LC-MS, Vaccine, Quality Control, Identity Testing, Regulation

Introduction

A Liquid Chromatography-Mass Spectrometry method (LC-MS) was developed, validated, and implemented in Quality Control (QC) for the identity testing of acellular pertussis combination vaccines. QC activities represent around 70% of product cycle time and are often considered as bottle necks. A single LC-MS method is able to replace several antibody based identity tests to enhance quality, increase efficiency, and reduce cost and cycle-time.

Methods

The vaccine samples are digested to generate peptides. The peptide mixtures are then separated using liquid chromatography and analyzed in real-time by the mass spectrometer. The mass spectrometer detects the antigen specific peptides or signature peptides using accurate mass and retention time.

Results

The LC-MS identity assay is able to not only detect the presence of each of the antigen specific peptides derived from Diphtheria Toxoid (D), Tetanus Toxoid (TT), Pertussis Toxoid (PT), Fimbriae (FIM), Pertactin (PRN), Filamentous Haemagglutinin (FHA), Inactivated Poliomyelitis Vaccine (IPV) but also the dosage of the vaccines for the identification and differentiation of the products with different formulations. The specificity of the method is validated by the positive detection of the antigen specific peptides in samples/products where the corresponding antigen is present and their lack of detection in samples/products where the corresponding antigen is absent. Robustness tests were performed to examine the potential sources of variation. This was done by selecting factors expected to affect method responses, and the effect of the factors were examined and statistically analyzed.

Conclusions

The identity testing for pertussis and related combination vaccines using LC-MS is specific and robust. The test method is validated and accepted by the Regulatory Agencies, such as FDA, Health Canada, and EU. Implementation of LC-MS ID test in QC replaces identity tests using traditional Western Blot, Dot Blot, Agar Immunodiffusion, and ELISAs, eliminates the use of more than 18 antibody reagents, and reduces the cost and cycle time.

Novel Aspect

Implementation of a LC-MS release test for vaccine identification in QC and obtaining regulatory approval represent breakthroughs in the vaccine industry.

WOr.05 - Food integrity, authenticity & traceability

Chairs: Ruud Peters, Michele Suman

WOr.05-1 Keynote: A STORY ON FOOD, FAITH AND FRAUD – DETECTION OF THE FAKES

Saskia van Ruth

Wageningen University & Research

Keywords: Authentication, criminals, food, fraud control

Food fraud is a form of criminal behaviour, no matter the definition of crime. Its consequences are devastating. The interaction between motivated offenders, and the opportunities presented by victims and lack of control measures favour occurrence of food fraud. Control measures help to counteract fraud opportunities and motivations. Analytical testing is one of those control measures. Traditional measurements have focused on the analysis of one or a few product characteristics. However, nowadays analytical techniques generating detailed analytical fingerprints are used to determine the identity of foods, and many techniques are available. Mass spectrometry methods are an important sub-group of these methods. Sensitive targeted methods will be presented, for instance for authentication of olive oil grades, as well as various fingerprint approaches. They include for instance broad anomaly methodology for authentication of spices, as well as for the characterization of organic products. Furthermore, MS-based breath analysis will demonstrate the importance of the origin of cocoa beans for a distinct identity of chocolates.

Conclusions

Food fraud is everywhere and its management requires analytical methods for detection of anomalies and confirmation of the authenticity of food products. Mass spectrometry offers a wide range of high value authentication applications.

WOr.05-2 DETECTING FOOD FRAUD: DETERMINATION OF A POTATO'S VARIETY AND ITS GEOGRAPHICAL ORIGIN IN ONE ANALYSIS VIA UPLC-IMS-QTOF

Christin Claassen ^(1,2) - **Juergen Kuballa** ⁽¹⁾ - **Sascha Rohn** ⁽²⁾

⁽¹⁾ GALAB Laboratories GmbH, Research & Development, Hamburg (Germany) - ⁽²⁾ Hamburg School of Food Science, University of Hamburg, Institute of Food Chemistry, Hamburg (Germany)

Keywords: Food fraud, metabolomics, lipid profiling, potato variety, geographical origin

Introduction

In recent years, consumers increasingly claim for regional potatoes of certain varieties. Suitably labelled products are considered to be of higher quality and prices what might result in increased rates of food fraud.

While regulators enhance consumer protection, existing analytical methods cannot ensure compliance of food properties and its labelling sufficiently. Therefore, a Metabolomics-based approach was developed to evaluate potatoes' authenticity.

Methods

About 200 organic and conventionally produced potato samples of known variety were collected from various German regions in 2016 and 2017.

Samples were cryomilled, lyophilized and extracted with a liquid-liquid-extraction protocol including chloroform, methanol and water. The analysis of nonpolar extracts was performed by an UPLC-IMS-QToF system. Obtained data sets were evaluated via multivariate data analysis using commercial and open source software.

Results

PCA and OPLS-DA analysis were carried out to differentiate between the varieties Belana, Gunda, Wega and Queen Anne. Additionally, samples were distinguished according to their growing regions. The structural elucidation of biomarkers was performed considering the retention time, precise molecular mass, MS/MS fragmentation data, isotopic patterns and ion mobility spectra provided by the UPLC-IMS-QToF analysis platform.

A set of more than 60 lipid biomarkers was identified that allows a reliable discrimination of the four potato varieties on the one hand and their cultivation region on the other. Characteristic substances, predominantly tri- and diglycerides, were found to differentiate potatoes originating from Thuringia, North Rhine-Westphalia, Bavaria, Schleswig-Holstein and other German production areas.

These generic markers withstand impacts like different farming methods and multiple harvesting periods. Furthermore, the reliability of the model was tested by analyzing mixtures of varieties and geographical origins.

Conclusions

The results demonstrate that this metabolomics-based approach in conjunction with a comprehensive biomarker database has potential as a future screening tool for the investigation of labelled potato tubers regarding their variety and their German cultivation region to prevent food fraud. This tool may be used to preserve consumer interests and support regulators' efforts to keep pace with fraudsters.

Novel Aspect

This approach determines a potato's variety and its cultivation region within one analysis whereas methods like electrophoresis and isotope analysis require multiple tests.

WO.r.05-3 NOSESPACE ANALYSIS BY PTR-TOF-MS FOR THE EXTRA VIRGIN OLIVE OIL CHARACTERIZATION

Julia Khomenko ^(1,2) - ***Eugenio Aprea*** ⁽¹⁾ - ***Flavia Gasperi*** ⁽¹⁾

⁽¹⁾ *Fondazione Edmund Mach, Research and Innovation Center, Department of Food Quality and Nutrition – via E. Mach 1, 38010 San Michele all'Adige (TN, Italy)* - ⁽²⁾ *Institute for Ion Physics and Applied Physics, University of Innsbruck, Technikerstr. 25, Innsbruck (Austria)*

Keywords: *direct injection, volatile organic compounds, inter-individual differences*

Introduction

Extravirgin olive (EVO) oil aroma and taste influence consumer preferences. Volatile organic compounds (VOCs) release from food within oral cavity is a complex process. Nose-space (NS) analysis should unravel the retro-nasal olfaction processes happened during olive oil tasting. Proton Transfer Reaction-Mass Spectrometry with Time-of-Flight analyser (PTR-ToF-MS) was already successfully applied for NS analysis of apples [2], cereal bars [3], coffee [4,5].

Methods

In this study 2 Italian EVO oils were used: commercial vs home made (Lago di Garda, TN). Pure oil, oil mixed with tomato sauce, oil on bread prepared for each EVO oil were used for headspace (HS) and NS measurements by PTR-ToF-MS 8000. The samples were presented to 8 assessors in duplicate. For NS sampling a glass nose-piece with silicone rubber tube was fitted into assessors' nostrils and connected to PTR-ToF-MS. The method for HS analysis was described elsewhere [6].

Results

26 out of 319 mass peaks of NS dataset displayed signal related to olive oil, tomato or bread aroma. Peaks related to human metabolism were eliminated. The individual profiles differed in

peak intensity and duration. However it was possible to distinguish the two types of EVO oils. HS data confirmed diversity of the oils. Home-made EVO oil in all products was characterized by significantly higher concentrations of mass peak 99.081 tentatively identified (t.i.) as (E)-2-hexenal associated to a green-leaf and fruity flavor [1]. Commercial EVO oil was richer in t.i. pentanal, acetic acid, butanal, and several esters which are reported as oxidation products of olive oil [7]. The samples of olive oil on bread showed the best performance (higher signals intensities) for all panellists due to an increased release surface area, prolonged presence in mouth, increased interaction with saliva and mucosal surfaces. Adding tomato sauce to olive oil decreased the olive oil aroma release in NS and contaminated it with tomato aroma.

Conclusions

NS analysis performed by PTR-ToF-MS discriminated different EVO olive oils. HS and NS data were consistent. The next steps are the optimization of sample preparations, tasting protocol, and the improvement of PTR-ToF-MS sensitivity by applying an ion funnel.

The optimized methods will be applied to study different EVO olive oils within the project "VIOLIN" aiming to the valorisation and promotion of Italian EVO oils.

Novel Aspect

For the first time NS analysis was applied for olive oil profiling taking into consideration individual perception and the combination with different food matrices.

References

1. C. M. Kalua, M. S. Allen, D. R. Bedgood Jr, A. G. Bishop, P. D. Prenzler, K. Robards; *Food Chemistry*, 100 (2007), pp.273-286
2. V. Ting, C. Soukoulis, P. Silcock, L. Cappellin, A. Romano, E. Aprea, P. Bremer, T. Märk, F. Gasperi, F. Biasioli; *Journal of Food Science*, 77(11) (2012), pp 1226-1233
3. S. Heenan, C. Soukoulis, P. Silcock, A. Fabris, E. Aprea, L. Cappellin, T. Märk, F. Gasperi, F. Biasioli; *Food Chemistry* 131 (2012), pp 477-484.
4. A. Romano, L. Cappellin, V. Ting, E. Aprea, L. Navarini, F. Gasperi, F. Biasioli; *International Journal of Mass Spectrometry*, 365-366 (2014), pp 20-27.
5. M. Charles, A. Romano, S. Yener, M. Barnabà, L. Navarini, T.T. Märk, F. Biasioli, F. Gasperi; *Food Research International*, 69(2015), pp 9-20.
6. S. Yener, A. Romano, L. Cappellin, T. Märk, J. Sánchez del Pulgar, F. Gasperi, L. Navarinic, F. Biasioli; *Journal of Mass Spectrometry*, 49 (2014), pp 929–935
7. M. Morales & R. Przybylski; *Olive Oil Oxidation in Handbook of olive oil. Analysis and properties*, Springer Science+Business Media New York 2013

Acknowledgment

The project was funded by Cariplo Foundation within the "Agroalimentare e Ricerca" (AGER) program. Project AGER2-Rif.2016-0169, "Valorizzazione dei prodotti italiani derivanti dall'oliva attraverso Tecniche Analitiche Innovative" - "VIOLIN".

WOr.05-4 **ULTRA-FAST IDENTIFICATION OF FISH SPECIES BY PAPERSPRAY COUPLED TO A BENCHTOP ORBITRAP MASS SPECTROMETER**

Daniel Lopez-Ferrer - Greg Foster - Mari Prieto

Thermo Fisher Scientific, Vertical Marketing, San Jose

Keywords: Food authenticity, PaperSpray, Orbitrap technology, Automations

Introduction

The largest seafood fraud investigation in the world to date found that 33% of seafood samples

are mislabeled in the U.S. Seafood fraud hurts consumers and also raises serious preservation concerns for our oceans. DNA-based analytical techniques are commonly used to determine the species that belongs to a given sample. However, the long analysis time complicates its applicability to be enforced. Here, we present a protein based method that identifies the most common commercially available salmon species in minutes using a simple protein extraction protocol coupled to PaperSpray® and top-down proteomics. The high mass accuracy and resolution of the benchtop quadrupole orbitrap MS allows not only the classification of proteins with high protein sequence homology, but also species dependent amino acid substitutions. The use of PaperSpray, a direct ionization technique, enabled us to simplify the mass spectrometric analysis and subsequent identification of the target proteins.

Methods

Seven reference samples from commercial species and 50 commercial foodstuffs were included in the work. Protein extraction was carried out by homogenizing 1 g of muscle. Water soluble proteins were centrifuged, the supernatant heated at 70 °C for 5 min and centrifuged again. Soluble proteins then were cleaned using automated microSPE devices. The analytes were then spotted onto the PaperSpray cassette and placed into the ProSolia Velox™ 360 PaperSpray® system coupled to a Thermo Scientific Q Exactive HF mass spectrometer. PaperSpray technology relies on cellulose paper and ions are generated directly from paper when an applied high voltage induces electrospray from the sharp tip of the wet paper. Proteins migrate through the paper and are directly ionized into the mass spectrometer. Extraction solvent used was 95/5/0.01 methanol/100 mM ammonium acetate/formic acid. Full MS spectra of intact thermostable proteins were analyzed using Thermo Scientific Biopharma Finder 2.0 software for molecular mass determination and the protein identities were confirmed by MSMS.

Results

The main goal of this work was to develop a simple and fast strategy that allows the authentication of fish species and be widely applied to enforce the guidelines from the U.S. Food and Drug Administration. The work was divided in two phases; primarily a discovery phase to develop the analytical workflow and to choose the best protein biomarker from standardized samples. This phase was followed by a validation phase to determine whether a commercial sample is mislabeled or not. During the discovery phase, a two-step workflow was created. First, the fraction of thermo-stable proteins from the tissue was purified because most of the commercial foodstuff has already been cooked. The second step, consisted of an ultrafast top-down proteomic approach of the purified proteins. Samples were cleaned up using reverse phase micro-spe cartridges. Proteins were eluted using 50% acetonitrile enabling easy spotting into the paper spray cassette while significantly reducing the drying time before being electro-sprayed into the mass spec. The overall sample preparation workflows allows for automated sample preparation and clean-up of the samples before been introduced in the MS system without impacting the throughput. A 1 min acquisition was performed and the proteins were then analyzed doing a top-down msx in the HCD cell of a benchtop quadrupole orbitrap mass spectrometer. For all the samples, an 11 kDa protein was the most abundant protein and identified as a parvalbumin (PRVB). This protein is highly conserved between different species. However, slightly mass differences in the mass of the PRVBs corresponding to amino acid substitutions allow for unambiguous systematic discrimination of the different species. Further validation on commercial samples demonstrated that over 10% of the samples were mislabeled. We believe that this workflow coupled to the PaperSpray technology platform could represent the ultimate payoff to reliable and widely used top down proteomics to trace and authenticate food.

Conclusions

In this work, we coupled PaperSpray ionization with high resolution accurate mass (HRAM)

Orbitrap technology for authentication of fish products. PaperSpray demonstrated that qualitative analysis of proteins is possible without time consuming sample preparation and advance chromatography techniques. This work also demonstrates the potential use of PaperSpray coupled to mass spectrometry as an emerging approach to protect customer rights and preventing fraudulent practices. And that can be useful to food industries that are seeking the opportunity to assure their food products labeling.

Novel Aspect

Fish traceability in minutes using a high throughput PaperSpray-Orbitrap MS platform and top-down proteomics.

WOr.05-5 CHARACTERIZATION AND CLASSIFICATION OF SPANISH PAPRIKA WITH PROTECTED DESIGNATION OF ORIGIN BY UHPLC-HRMS POLYPHENOLIC FINGERPRINTING AND CHEMOMETRICS

Sergio Barbosa - Guillem Campmajó - Oscar Núñez - Javier Saurina - Lluís Puignou

University of Barcelona, Department of Chemical Engineering and Analytical Chemistry, Barcelona

Keywords: Polyphenols; Spanish Paprika; Food Authentication; UHPLC; High resolution mass spectrometry

Introduction

Paprika is a red powder seasoning with a characteristic flavor obtained from red peppers. Spanish paprikas from La Vera and Murcia have protected designation of origin (PDO), which increases product value but also their production costs, making necessary to properly authenticate them. Polyphenolic contents in these products may be employed as descriptors to achieve their authentication [1].

Methods

A UHPLC-HRMS (Orbitrap) method was employed to obtain paprika polyphenolic fingerprints. HRMS data was processed with TraceFinder v3.3 software by employing a customized target accurate mass database of 54 polyphenols. Polyphenolic data was considered as a source of potential descriptors for the classification and characterization of paprika by exploratory principal component analysis (PCA).

Results

111 samples belonging to La Vera and Murcia Spanish PDOs (including sweet, bittersweet and spicy varieties), as well as 15 samples obtained from the Check Republic were analyzed. The plot of scores obtained after PCA using UHPLC-HRMS polyphenolic profiles revealed patterns that were correlated to sample characteristics such as the production country (Spain, Check Republic) and regions, allowing a clear differentiation between paprika PDOs and varieties. Based on the PCA plot of loadings, the most remarkable polyphenols and phenolic acids on each type of paprika PDO and variety were identified, making possible to differentiate them according to their different content and distribution of polyphenols. These polyphenols can be studied as possible future biomarkers to guarantee sample authentication and for the prevention of frauds.

Conclusions

The obtained results demonstrate that UHPLC-HRMS polyphenolic profiles obtained by a simple screening with a customized accurate mass database of polyphenols can be employed to achieve the characterization and classification of paprika samples according to their production region to guarantee product authentication. Potential polyphenolic bioactive compounds for paprika varieties were identified.

Novel Aspect

New analytical approach for the characterization and classification of Spanish paprika based on HRMS polyphenolic fingerprinting.

References

1. Lucci P., Saurina J., Núñez O., TrAC Trends Anal. Chem., 88, 1 (2017).

WOr.06 - MS tools in doping & toxicology

Chairs: Mario Thevis, Gianluca Bartolucci

WOr.06-1 **Keynote:** EVOLUTION AND APPLICATION OF MASS SPECTROMETRY –BASED ANALYTICAL METHODS IN ANTI-DOPING ANALYSES

Tiia Kuuranne

Swiss Laboratory for Doping Analyses, University Center of Legal Medicine Geneva and Lausanne, Centre Hospitalier Universitaire Vaudois and University of Lausanne, Lausanne, Switzerland

Keywords: anti-doping, WADA accredited laboratory, prohibited substances, sport, mass spectrometry

Introduction

Accredited anti-doping laboratories operate within strict regulations of WADA, such as the International Standard for Laboratories (ISL) and technical documents, which aim at harmonizing performance and reporting of routine analyses. Due to the sensitivity and selectivity required from confirmation analyses, mass spectrometry (MS) coupled to either gas (GC) or liquid chromatography (LC) is the analytical technique of choice for identification of the analytes.

Methods

Whereas the first routine applications were based on the GC/MS techniques, implementation of LC/MS-methodologies during the last 15-30 years dramatically extended the menu of prohibited substances amenable for MS-identification. Both platforms are, however, required for adequate coverage of analyses, which have mainly qualitative but also some quantitative targets, and benefit greatly of technical advances of sensitivity, resolution and data processing.

Results

In the evolution of the GC/MS-techniques, the cornerstones include the implementation of enhanced trimethylsilyl (TMS) reagent for keto-group, which had a significant impact on the robust analysis of not only exogenous, but also endogenous anabolic androgenic steroids (AAS) with keto-structure. Despite of roots originating to 1980's, this technique is still actively in use in steroid profiling for athlete biological passport (ABP) purposes. Regarding the endogenous AAS, such as testosterone, another major advance has been the development of an analytical approach where GC separation, combustion to CO₂, and isotope ratio mass spectrometry (GC/C/IRMS) was combined for the determination of doping use, i.e. exogenous origin of AAS. Along the LC/MS-based methods, especially those with soft ionization methods such as electrospray, anti-doping analyses have become more sophisticated for prohibited substances with low volatility (corticosteroids), instable derivatives (diuretics), and higher molecular mass (conjugates, peptides).

Conclusions

The spectrum of target compounds in anti-doping analyses is wide, dynamic and strongly

dependent on the pharmaceutical drug development. In order to keep up with the analytical capabilities for emerging performance enhancing substances, the laboratories should maintain the flexibility of analytical pipelines but also follow the development of mass spectrometers and their performance to update the processes for improved sensitivity, selectivity and robustness.

Novel Aspect

Samples can be stored up to ten years. With re-analyses, potential of non-restrictive data acquisition and digital matrix are proposed to increase deterrence effect of doping control.

WOr.06-2 LC-HRMS/MS VS. LC-MRM: SCREENING OF COVALENT DRUG-PROTEIN ADDUCTS

Timon Geib - Lekha Sleno

Université du Québec à Montréal, Chemistry Department, Montréal

Keywords: clozapine, acetaminophen, covalent protein binding, multiple reaction monitoring, high resolution tandem mass spectrometry

Introduction

Hepatotoxicity of prescription and over-the-counter medication is resulting in major drug-induced adverse effects. Acetaminophen (APAP) and clozapine (CLZ) are known to form reactive metabolites, which can covalently bind to cysteine residues in hepatic proteins. We have investigated the protein binding of these two compounds to different glutathione S-transferase (GST) isoforms.

Methods

Our developed method uses *in vitro* activation of APAP and CLZ to NAPQI and CLox, respectively, with either liver microsomes or CYP3A4 Supersomes, while adding purified GSTs to the incubation. Further sample preparation employed tryptic or peptic digestion, coupled to LC-MS/MS analysis. Results between SPE, immunoprecipitation and offline LC fractionation have been compared. We have developed two LC-MS/MS based strategies for screening for modified peptides.

Results

Method optimization was carefully performed on different cytosolic *Homo sapiens* GSTs. LC-MS/MS of information dependent acquisition (IDA) files were searched with UniProtKB/Swiss-Prot protein database. The search algorithm was altered to find NAPQI- or CLZox-cysteine adducts. Multiple reaction monitoring (MRM) was also used as a targeted screening technique. Optimal conditions for incubations were determined (3 h at 37 °C), followed by either tryptic or peptic digestion at 37 °C for 4 h or 1 h, respectively.

The first LC-MS/MS screening technique used a traditional shotgun proteomics workflow (IDA) on a QqTOF with subsequent database searching. Then, a targeted MRM method for all cysteine-containing peptides of the individual GST isoforms was tested. The GSTM2 and P1 were found to be a target of NAPQI *in vitro*, with more sites identified using the MRM based screening method. Cys174 in M2 and Cys48 in P1 were found using IDA and MRM; while additional sites Cys115 in M2, and Cys15 and Cys170 in P1 were confirmed by MRM only.

Conclusions

Currently, CLZox-protein binding, as well as, IP purification are being optimized to provide highly purified NAPQI- and CLZox-peptide extracts for fast LC-MS/MS analyses of complex samples, including exposed liver microsomes, S9 fractions and homogenates. Comparison of the performance of IDA and MRM analyses showed that MRM represents a complementary

technique to untargeted IDA, which can be easily applied in future drug-protein binding investigations.

Novel Aspect

High sensitivity screening and identification of drug-related protein binding to glutathione S-transferases using LC-HRMS/MS and LC-MRM

WOr.06-2 COUPLING MASS SPECTROMETRY AND CHEMOMETRICS FOR STRUCTURAL CHARACTERIZATION OF ANDROGENIC ANABOLIC STEROIDS

Patrizia Leogrande⁽¹⁾ - **Xavier de la Torre**⁽¹⁾ - **Francesco Botrè**^(1,2) - **Daniel Jardines**⁽¹⁾ - **Federico Marini**⁽³⁾

⁽¹⁾ Laboratorio Antidoping, Federazione Medico Sportiva Italiana, Largo Giulio Onesti 1, 00197, Rome, Italy -

⁽²⁾ Dipartimento di Medicina Sperimentale, "Sapienza" Università di Roma, Viale Regina Elena 324, 00161 Rome, Italy - ⁽³⁾ Dipartimento di Chimica, "Sapienza" Università di Roma, Viale Regina Elena 324, 00161 Rome, Italy

Keywords: Anabolic Androgenic Steroids, Doping Control, Mass Spectrometry, Multivariate Analysis, Predictive Models

Introduction

The Androgenic Anabolic Steroids (AAS) are one of the main and most used classes of doping substances, included in the S1.1 section of the prohibited list of the World Anti-Doping Agency (WADA)[1].

An increasingly widespread phenomenon is that of designer steroids, molecules designed and synthesized by modifying the structure of known AASs, in the aim of maintaining the same biopharmacological activity but resulting undetectable at the anti-doping test.

Methods

Predictive models were constructed based on the MS characterization of 133 AAS compounds, to identify new steroids and obtain information on their molecular structure. Starting from Principal Components Analysis different chemometric methods were applied, such as classification and clustering methods, outlining a spectral and structural characterization for each subclass of AAS

Results

Mass spectrometric data obtained by GC-qTOF experiments performed on the TMS-derivatives allowed to group the AAS in 5 subclasses. PCA showed the high similarity between the class of 4-ene-3-keto steroids and the steroids substituted on the A ring.

A discriminant classification step by PLS-DA was performed to build a predictive model to separate the subclasses and to identify the most significant variables, like the unsaturation in position 17 on D ring, characteristic of the 3-hydroxy steroids category.

The use of a modeling classification method, which operates by building individual models for the different categories, allowed to define the characteristics of uniqueness of the investigated steroid subclasses. This approach resulted particularly effective to identify the salient features of the class of 5-ene steroids.

Clustering methods were also applied, with results consistent with most of the previously obtained outcomes.

All the information has been summarized by identifying class specific common fragments and spectral trends.

Conclusions

The results of this study suggest that it is possible to trace, as quickly as possible and with the

highest degree of certainty, any novel doping substances – and especially “designer AAS” – illicitly synthesized in such a way as to be invisible to the current tests.

The newly developed methods to reduce the time for the identification of unknown substances by fitting/unfitting them in one of the specific subclasses.

Novel Aspect

The predictive nature of the chemometric models was based on the characterization of 133 AAS, and considering the contribution of more than 30 variables.

References

WADA. The World Anti-Doping Code International Standard. Prohibited List, January 2018. Available at: https://www.wada-ama.org/sites/default/files/prohibited_list_2018_en.pdf.

WOr.06-4 DETECTION OF CHEMICAL WARFARE AGENTS-RELATED PHENYLARSENIC COMPOUNDS IN MARINE BIOTA AND IDENTIFICATION OF THEIR METABOLITES IN IN VITRO STUDIES BY ORBITRAP HRMS

Hanna Niemikoski⁽¹⁾ - ***Daniel Koske***⁽²⁾ - ***Paula Vanninen***⁽¹⁾

⁽¹⁾ Finnish Institute for Verification of the Chemical Weapons Convention, University of Helsinki, Helsinki - ⁽²⁾ Thünen Institute, Institute of Fisheries Ecology, Hamburg

Keywords: Chemical warfare agent, arsenic, fish, metabolism, high resolution mass spectrometry

Introduction

Phenylarsenic compounds, Clark I, Clark II and Adamsite, were designed as riot control agents during World War I [1]. Nowadays, these chemicals are classified as chemical warfare agents (CWAs) by Chemical Weapons Convention. These agents were produced in a large scale during World War I and II. After World War II, their disposal began by sea-dumping. Hundreds of thousands of tons of chemical ammunitions (mainly sulfur mustard and phenylarsenic compounds) were dumped in the Baltic Sea and Skagerrak area between Norway and Denmark [2]. Several investigations have shown that sediment samples collected near these dumping areas are contaminated with CWAs as a result of leaking containers [3]. These leaking toxic chemicals are posing a threat to marine environment. In aqueous environment phenylarsenic compounds degrade forming hydrolysis and oxidation products. Previously traces of oxidation form of Clark I and/or II have been detected from lobster and flatfish samples collected from dumpsite area near Swedish coastline [4]. In order to study uptake of arsenic containing CWAs in fish in larger scale, fish samples were collected from different dumpsites and analysed. One goal of the ongoing EU-project DAIMON (Decision Aid for Marine Munitions) (2016-2019) is assessment of munitions' impact on marine biota. One aim is to develop specific and sensitive methods to assess the fate and impact of toxic CWAs on marine biota [5]. Quantitative analysis of these chemicals and their degradation products are needed to prove presence of these products in aquatic biota in order to support risk assessment for possible accumulation in food chain. There are hardly any studies about transference of arsenic based CWAs into aquatic biota. Due to the lack of information how these agents behave in living organism it's impossible to predict interactions or transformations of these chemicals. In order to study transformations of arsenic containing CWAs in fish we performed metabolism studies with fish liver homogenates. Information achieved from metabolism studies is essential when developing novel analysis methods for arsenic containing CWAs from fish samples.

Methods

Fish samples obtained from three different dumpsite areas were analysed by liquid chromatography-tandem mass spectrometry (LC-MS/MS) in the selective reaction monitoring

(SRM) mode. Qualitative identification criteria based on EU guidelines were used to reliable identification of detected CWA-related compounds. In order to study transformations of phenylarsenic compounds in fish *in vitro* we incubated five CWA-related arsenic compounds in cod liver homogenates (S9 fraction). Protein concentration of liver homogenates was determined by using the Lowry method. Samples from metabolism studies were analysed by liquid chromatography-Orbitrap high resolution mass spectrometry.

Results

Fish samples collected from three different dumpsites contained CWA-related arsenic compounds at parts per billion (ppb, ng/g) concentration levels. These arsenic containing chemicals were analysed as their pentavalent forms by LC-MS/MS method. Metabolism studies with cod liver homogenates proved that phenylarsenic compounds as trivalent intact form as well as oxidized pentavalent form formed glutathione (GSH) conjugates. Structures of these metabolites were determined by high resolution mass spectrometry. Also reference chemicals of detected glutathione metabolites were synthesized and their structures were confirmed additionally by nuclear resonance spectroscopy. MS spectra of detected metabolites from *in vitro* metabolites studies were found to match spectra of synthesized reference chemicals.

Conclusions

There are hardly any studies about transference of arsenic based CWAs into aquatic biota. Previously these arsenic-containing chemicals have been detected from marine biota samples as their pentavalent oxidation state. Metabolism studies with fish liver enzymes proved that these oxidized arsenic compounds were reduced by glutathione and formed glutathione conjugates. There are no data available about toxicities of these trivalent phenylarsenic-glutathione metabolites neither how these metabolites behave in fish. It's possible that these glutathione metabolites are excreted into bile or urine. Based on this, more studies are needed to evaluate total arsenic-related CWAs concentrations in aquatic biota.

Novel Aspect

CWA-related arsenic compounds have been detected in different marine species for the first time. In order to get more specific information about transformations of these arsenic compounds in fish, metabolism studies with fish liver homogenates were applied. Phenylarsenic compounds formed glutathione conjugates which structures were confirmed by Orbitrap high resolution mass spectrometry. None of these metabolites have ever been detected from marine biota samples and this is the first time when metabolism of phenylarsenic compounds have been studied with fish liver homogenates *in vitro*.

References

1. Szinicz L., Toxicology, 214,67–181(2005).
2. Glasby G., Sci.TotalEnvironm.,206, 267–273 (1997).
3. Bełdowski J., Klusek Z., Szubska M., Turja R., Bulczak A.I., Rak D., Brenner M., Lang T., Kotwicki L.,Grzelak K., Jakacki J., Fricke N., Östin A., Olsson U., Fabisiak J., Garnaga G., RattfeltNyholm J., Majewski P., Broeg K., Söderström M., Vanninen P., Popiel S., Nawala J., Lehtonen K., Berglind R., Schmidt B., *Deep-Sea Res Part II*, 128, 85–95 (2016).
4. Niemikoski H., Söderström M. and Vanninen P., *Anal. Chem.* 89, 11129–11134 (2017).
5. <http://www.daimonproject.com>

WOR.06-5 A FIT-FOR-PURPOSE UHPLC-MS/MS METHOD FOR ROUTINE SCREENING OF SARM RESIDUES IN BOVINE AND EQUINE SERUM

***Emiliano Ventura*⁽¹⁾ - *Anna Gadaj*⁽¹⁾ - *Tom Buckley*⁽²⁾ - *Mark H. Mooney*⁽¹⁾**

⁽¹⁾ Institute for Global Food Security, Queen's University Belfast, School of Biological Sciences, Belfast - ⁽²⁾ Irish Equine Centre, Irish Diagnostic Laboratory Services Ltd., Johnstown, Co. Kildare

Keywords: SARMs, UHPLC-MS/MS, serum, doping control, food safety

Introduction

Selective Androgen Receptor Modulators (SARMs) are an attractive alternative to anabolic-androgenic steroids despite the fact they are not yet approved as pharmaceutical drugs [1]. Consequently, there is much evidence of their availability through black and grey market sources. In this study, a semi-quantitative UHPLC-MS/MS method was developed to monitor the abuse of 15 compounds belonging to eight different SARM families, in bovine and equine animals, respectively.

Methods

Serum samples underwent protein crash (0.5 mM NH₄OH in MeCN), induced phase separation (NaCl) and defatting (*n*-hexane) before concentration under nitrogen. Extracted SARM residues were subsequently analysed by UHPLC-MS/MS operating in ESI+/ESI-. Chromatographic separation was performed using a Luna Omega Polar C18 column at 45°C, and a binary gradient elution of 14 minutes employing H₂O and MeOH, both containing 0.1% (v/v) CH₃COOH as mobile phases at a flow rate of 0.4 mL/min.

Results

Validation of the semi-quantitative assay was carried out according to the Community Reference Laboratories Residues (CRLs) 20/1/2010 guidelines [2]. The following performance studies were carried out during the validation process: specificity, selectivity, detection capability (CC_β), absolute recovery as well as applicability, ruggedness and matrix effects. The detection capability was calculated by assessing threshold value (T) and cut-off factor (F_m) and CC_β values were determined at the following concentrations: 0.5 ng mL⁻¹ (CDS-025139), 1 ng mL⁻¹ (AC-252636, PF-06260414), 2 ng mL⁻¹ (bicalutamide, GLPG-0492, LGD-2226, S-1, S-6, S-22, S-23) and 5 ng mL⁻¹ (BMS-564929, LGD-4033, , RAD140, S-4, S-9), respectively. Moreover, ongoing work aims to identify *in vivo* metabolites of selected SARMs compounds following animal administration studies, which could serve as additional marker residues to be included within the current test method.

Conclusions

SARMs are widely recognised as emerging drugs of abuse in animal sport, and also as potent candidates for misuse in food producing animals. The developed method is currently being employed as a rapid, simple, environmental friendly and cost-effective tool for the screening of serum samples to evaluate the abuse of SARM drugs in stock farming and to monitor for doping practices in racing horses, ensuring customer safety and fair play in animal sport competitions.

Novel Aspect

The developed high-throughput assay allows, through a relatively short run and a small test sample volume (400 μL) of serum, the analysis of a total of 50 samples in a single day.

References

1. Thevis M., Schänzer W., *Molecular and Cellular Endocrinology*, 464, 34-45 (2018).
2. Community Reference Laboratory, *Guidelines for the validation of screening methods for residues of veterinary medicines* (2010).

Acknowledgements

The research was supported by funding from the European Union Horizon 2020 Research and Innovation Action under Grant Agreement No. 642380.

WOr.07 - MS imaging: applications

Chairs: Martina Marchetti-Deschmann, Fulvio Magni

WOr.07-1 **Keynote: COMBINING MASS SPECTROMETRY IMAGING & MICROPROTEOMICS TO INVESTIGATE INTRATUMOR HETEROGENEITY**

Liam McDonnell

Leiden University Medical Center

Keywords: MALDI, imaging, microproteomics, laser capture microdissection, intratumor het

Introduction

MSI simultaneously records the distributions of hundreds of molecules directly from tissue. When combined with multivariate analysis it can reveal regions of tissue with distinct molecular signatures, including tumor subpopulations.¹ The small pixel sizes and direct tissue analysis of MSI limits the achievable depth of proteome coverage. Here, we demonstrate how MSI can be combined with high sensitivity microproteomics, even of the same tissue section.

Methods

Protein MALDI MSI data of patient tissues was combined with histopathological analysis to isolate the MSI data from tumor regions. Five independent multivariate tests were used to identify tumor subpopulations, and subpopulations associated with phenotype identified. Regions of tissue were isolated by laser capture microdissection and analyzed by LC-MS/MS using an Orbitrap Fusion.

Results

MALDI MSI has been applied to patient tissues of gastric cancer, breast cancer, and soft tissue sarcomas. Following histopathological annotation and automated feature identification and extraction, the resulting datacube of MS images of only tumor regions, was subject to five independent multivariate methods to identify tumor subpopulations characterized by distinct molecular signatures. Associations with clinical data identified those subpopulations linked with patient phenotype (survival and metastasis).

In order to dig deeper into proteome changes in small, localized regions of tissue we have developed an SP3 based microproteomics approach utilizing 10-plex TMT labels and high-pH fractionation, to combine high relative quantitation precision with in-depth proteome profiling. We demonstrate the performance of this approach and demonstrate how it may be seamlessly combined with MALDI MSI of the same tissue section.

Conclusions

MALDI MSI combined with microproteomics can be used to identify and investigate phenotypic tumor subpopulations.

Novel Aspect

Intratumor heterogeneity; seamless integration of MALDI MSI and LC-MS/MS, of the same tissue section.

References

1. Balluff B., Frese C.K., Maier S.K., Schöne C., Kuster B., Schmitt M., Aubele M., Höfler H., Deelder A.M., Heck A.J.R., Hogendoorn P.C.W., Morreau J., Altelaar A.F.M., Walch A., McDonnell L.A. *J. Pathol.* 235, 3-13 (2015).
2. deGraaf E.L., Pellegrini D., McDonnell L.A. *J. Proteome Res.* 15, 4722–4730 (2016).

WOR.07-2 MULTISITE MULTIMODAL MASS SPECTROMETRY IMAGING OF GENETICALLY ENGINEERED MOUSE MODELS CONTAINING MUTATIONS LINKED TO BOWEL CANCER

Chelsea J. Nikula⁽¹⁾ - **Rory T. Steven**⁽¹⁾ - **Efstathios Elia**⁽¹⁾ - **Andrew D. Campbell**⁽²⁾ - **Spencer A. Thomas**⁽⁷⁾ - **Alex Dexter**⁽¹⁾ - **Teresa Murta**⁽¹⁾ - **Stephanie Ling**⁽³⁾ - **Arafath K. Najumudeen**⁽²⁾ - **Adam Taylor**⁽¹⁾ - **Jean-Luc Vorng**⁽¹⁾ - **Gregory Hamm**⁽³⁾ - **Rasmus Havelund**⁽¹⁾ - **Renata Filipe-Soares**⁽⁴⁾ - **David Gay**⁽²⁾ - **Bin Yan**⁽¹⁾ - **Ala Al-Afeef**⁽¹⁾ - **Tingting Fu**⁽¹⁾ - **Robin Philip**⁽¹⁾ - **Kenneth N. Robinson**⁽¹⁾ - **Simon Barry**⁽³⁾ - **Kevin Brindle**⁽⁵⁾ - **Ian S. Gilmore**⁽⁷⁾ - **Richard J. A. Goodwin**⁽³⁾ - **John F. Marshall**⁽⁶⁾ - **George Poulgiannis**⁽⁷⁾ - **Zoltán Takáts**⁽⁴⁾ - **Mariia O. Yuneva**⁽⁸⁾ - **Owen J. Sansom**⁽²⁾ - **Josephine Bunch**⁽¹⁾

⁽¹⁾ National Physical Laboratory, NiCE-MSI, Teddington - ⁽²⁾ Cancer Research UK Beatson Institute, Department of Invasion and Metastasis, Glasgow - ⁽³⁾ AstraZeneca, iMED, Cambridge - ⁽⁴⁾ Imperial College London, Department of Surgery and Cancer, London - ⁽⁵⁾ Cancer Research UK Cambridge Institute, University of Cambridge, Cambridge - ⁽⁶⁾ Barts Cancer Institute, Queen Mary University of London, London - ⁽⁷⁾ Institute of Cancer Research, Division of Cancer Biology, London - ⁽⁸⁾ The Francis Crick Institute, Department of Tumour Biology, London

Keywords: mass spectrometry imaging, bowel cancer

Introduction

Bowel cancer is the second leading cause of cancer mortality in the United Kingdom. Mutations in the APC and KRAS genes are found in these tumours at rates of over 70% and 40%, respectively. In order to better understand the effects of these mutations on cellular metabolism, MALDI, DESI, and SIMS mass spectrometry imaging were used at multiple research sites to identify and map the localization of key metabolites in genetically engineered murine tissue.

Methods

Murine samples were prepared and sectioned onto glass or ITO coated slides and stored at -80°C until analysis. MALDI-MSI was carried out on a RapifleX ToF-ToF (Bruker Daltonics) or Synapt G2-Si Q-IMS-ToF (Waters). DESI-MSI was carried out on a Synapt or Xevo Q-TOF (Waters). SIMS analysis was carried out on an OrbiSIMS (IONTOF, Thermo Scientific). Data were converted to imzML by imzML converter and analysed in MatLab (2017a, MathWorks) via SpectralAnalysis.

Results

Colon and small intestine from sixteen animals were embedded into blocks of four containing four different genotypes each (wild type, KRAS-mutated, APC-mutated, and dual KRAS/APC-mutated). DESI and MALDI data were collected at multiple sites to discern variations in data seen between different instruments and to enhance method development to produce robust, reproducible practices. Preliminary analysis of these data indicates that the use of MSI followed by multiple computational analysis techniques including NMF, PCA, k-means clustering, and t-SNE neural networking allows for the segmentation of anatomical regions as well as grouping of phenotypic variation by genotype. Work is ongoing to utilize database matching and MS/MS studies to identify metabolites present in segmented data to correlate metabolic changes to genetic type. Lastly, further analysis of these datasets will allow guided analysis of additional samples via high resolution MSI to examine metabolite localization in a cellular context.

Conclusions

Preliminary analysis indicates that both anatomical and phenotypic segmentation can be generated using a combination of computational analysis techniques. Current work is underway to identify metabolites within these segmented regions and determine how the APC and KRAS mutations may be effecting metabolic processes. Future work can integrate MSI results with proteomics and transcriptomics to inform additional studies.

Novel Aspect

Multimodal, multisite analysis of bowel cancer murine models by MSI followed by multiple computational analysis methods to map metabolites using database matching.

WOr.07-3 CAN 3D MALDI-MS IMAGING BE THE MISSING PIECE FOR UNDERSTANDING THE TRAUMATIC BRAIN INJURY PUZZLE?

Khalil Mallah⁽¹⁾ - **Jusal Quanico**⁽¹⁾ - **Firas Kobeissy**⁽²⁾ - **Kazem Zibara**⁽³⁾ - **Michel Salzet**⁽¹⁾ - **Isabelle Fournier**⁽¹⁾

⁽¹⁾ INSERM, U1192 - Laboratoire Protéomique, Réponse Inflammatoire et Spectrométrie de Masse (PRISM), Université de Lille, Department of Biology, Lille - ⁽²⁾ American University of Beirut, Faculty of Medicine, Department of Biochemistry and Molecular Genetics, Beirut - ⁽³⁾ PRASE, ER045, Laboratory of stem cells, Lebanese University, Faculty of Sciences-I, Department of Biology, Beirut

Keywords: TBI, 3D MALDI-MSI, MALDI-LTQ-XL-Orbitrap

Introduction

Traumatic Brain Injury is mainly due to direct mechanical damage to the brain tissue leading to an axonal disruption and wide spread neural dysfunction. Monitoring lipid changes post injury could help in a possible therapeutic approach. In fact, lipid biomarkers have been studied in several brain injury models using mass spectrometry as demonstrated by Nielsen et al which identified NAPE as a biomarker of dead neurons in mice cerebral ischemia injury model.¹

Methods

Lipid 3-Dimensional MALDI mass spectrometry imaging was applied on injured brain harvested from rat subjected to controlled cortical impact model (CCI) injury model.² Images were conducted using high resolution MALDI-LTQ-XL Orbitrap. Data visualization and processing was performed on Scils lab MVS software along with 3D image reconstruction, principal component analysis, and unsupervised spatial segmentation.

Results

MALDI images were performed on serial sections from an injured rat brain at 3 days post injury and the 3D reconstruction was performed after performing unsupervised spatial segmentation on all sections followed by image alignment. The unsupervised spatial segmentation showed two unique clusters that co-localized within the injured cortical tissue and not expressed in any other regions of the brain. These clusters both departed from the non-injured cortical tissue on the opposite hemisphere of the brain in the PCA analysis, thus confirming the first finding of lesion specific clusters. In 3D manner, we visualized and assigned specific m/z values solely expressed in the injured tissue. At the same time, certain m/z ions showed a particular expression within the injury site and the ventricular system of the brain.

Conclusions

Our results can help in further monitoring and identifying lesion-specific m/z in a 3D manner to obtain a better understanding and visualization of molecular and cellular events occurring post-injury. We have also showed injury-elevated m/z in common with the ventricular system of the

brain thus suggesting a possible trafficking of molecules between the injured brain area and the spinal cord.

Novel Aspect

Applying 3D MALDI-MSI on traumatic brain injury for the first time. High spectral resolution allowed to differentiate between m/z ions with a difference of up to 0.1 mDa.

References

1. Nielsen, M. M. B.; Lambertsen, K. L.; Clausen, B. H.; Meyer, M.; Bhandari, D. R.; Larsen, S. T.; Poulsen, S. S.; Spengler, B. Mass Spectrometry Imaging of Biomarker Lipids for Phagocytosis and Signalling during Focal Cerebral Ischaemia. *Nat. Publ. Gr.* **2016**, No. May, 1–14.
2. Xiong, Y.; Mahmood, A.; Chopp, M. Animal Models of Traumatic Brain Injury. *Nat. Rev. Neurosci.* **2013**, *14*, 128–142.

WOr.07-5 SILVER-DOPED NANO-DESI MSI FOR PROSTAGLANDIN IMAGING

Kyle Duncan⁽¹⁾ - Ru Fang⁽¹⁾ - Jia Yuan⁽²⁾ - Rosalie Chu⁽³⁾ - Sudhansu K. Dey⁽²⁾ - Kristin Burnum-Johnson⁽³⁾ - Ingela Lanekoff⁽¹⁾

⁽¹⁾ Uppsala University, Dept of Chemistry-BMC, Uppsala - ⁽²⁾ Cincinnati Children's Hospital Medical Center, Division of Reproductive Sciences, Cincinnati - ⁽³⁾ Pacific Northwest National Laboratory, Biological Sciences Division, Richland

Keywords: nano-DESI MSI, quantification, prostaglandin, silver, uterine tissue

Introduction

Silver ions are known to form weak charge transfer complexes with alkenes,¹ and the addition of Ag⁺ to a double bond can therefore increase ionization efficiency and detectability of a compound. This work demonstrates nanospray desorption electrospray ionization mass spectrometry imaging (nano-DESI MSI) of endogenous prostaglandin species (PG) by utilizing the complexation between Ag⁺ and PG to enable the generation of unique PG ion images from thin tissue sections.

Methods

Nano-DESI MSI was performed by connecting a custom built source to a Q-ExactiveOrbitrap mass spectrometer. 10 ppm Ag⁺ was added to the solvent, which also contained acetonitrile and methanol (9:1), and deuterated PG species to facilitate quantification. Mouse uterine tissue sections of ~ 2 mm in diameter were imaged with higher spatial resolution using oversampling.² The generated quantitative ion images of PG had an approximate pixel size of 30 μm x 10 μm.

Results

The inclusion of silver into the solvent increased the sensitivity for [PG+Ag]⁺ ions by ~ 30 times compared to conventionally monitored [PG-H]⁻ ions. Additionally, it was found that the natural isotopic pattern from the two stable Ag⁺ isotopes was beneficial when corroborating the identity of analytes as silver ion adducts and that the exclusive binding of Ag⁺ to the alkene moiety enabled direct identification of double bond containing metabolite species.

Quantitative nano-DESI MSI revealed the localization of PG species to the luminal epithelium and glandular epithelium of mouse uterine tissue sections at day 4 of pregnancy prior to embryo attachment. Interestingly, the amount of both PGE₂ and PGF_{2α} were found to be three times as high in these cellular regions. The cellular regions of the glandular epithelium are only 100-300 μm in diameter, and were resolved in this study by use of oversampling. Three additional PG species were identified in the tissue with similar localizations displaying the great utility of silver-doped nano-DESI MSI.

Conclusions

The presented results successfully demonstrate quantitative nano-DESI MSI of [PG+Ag]⁺ species directly from biological tissue without additional sample preparation. The described approach will be extremely valuable for future studies investigating the role of PG species in biological systems.

Novel Aspect

First study detailing the specific localization and quantity of prostaglandin species in biological tissue sections.

References

1. Frenking, G.; Fröhlich, N. *Chem. Rev.* ,100, 717 (2000)
2. Duncan, K., Lanekoff I., *Anal. Chem.*, 90, 2451 (2018)

WOr.08 - Gas phase ion chemistry & thermochemistry

Chair: Pan Yanjiang

WOr.08-1 **Keynote:** PROGRESS TOWARDS ENZYME-LIKE SPECIFICITY IN PROTEIN FRAGMENTATION

Scott A. McLuckey - David J. Foreman

Department of Chemistry, Purdue University, West Lafayette, IN, USA 47907-2084

Keywords: top-down protein characterization, selective fragmentation, ion/ion reactions, aspartic acid effect, ornithine effect, dehydroalanine effect

Introduction

Strategies for top-down tandem mass spectrometry of protein ions have generally sought to maximize the number of structurally diagnostic product ions. We are interested in developing complementary approaches that tend to lead to a limited number of relatively abundant cleavages that arise from specific/predictable locations.

Methods

Experiments were conducted using a quadrupole/time-of-flight tandem mass spectrometer (Sciex 5600) that has been modified to conduct ion/ion reaction experiments. Protein precursor ions were subjected to ion/ion proton transfer reactions to bring them into a mass-to-charge region in which aspartic acid and proline cleavages are favored upon ion trap CID.

Results

Three strategies that provide enzyme-like specificity in protein ion fragmentation are described. These include taking advantage of the known tendencies for cleavages at aspartic acid and proline residues, as well as taking advantage of the ornithine and dehydroalanine effects. We describe an approach that concentrates precursor ion charge states derived from ESI under denaturing conditions into a single charge state via ion parking, ion isolation, proton transfer ion/ion reactions to generate precursor ions within a charge state range where aspartic acid and proline cleavages are most prominent, and use of a broad-band collisional activation technique to fragment all of the charge states simultaneously. This approach is demonstrated using four standard proteins commonly used to test/demonstrate top-down protein characterization approaches. The other two approaches rely on converting existing residues into either ornithine or dehydroalanine. These residues also give rise to highly favored cleavages.

Conclusions

The ultimate goal is to develop a suite of strategies that provides the analyst with options for multiple selective approaches for deconstructing mass-selected protein ions as complementary approaches to the standard means for maximizing cleavage in non-specific fashion. The present work demonstrates approaches that are selective for aspartic acid residues, proline residues, arginine residues, and several other residues that can be used to generate dehydroalanine.

Novel Aspect

Novel approaches are described that favor selective cleavages at a variety of residue types in protein ions.

WOr.08-2 DISSOCIATION OF HYDROGEN-ABUNDANT PEPTIDE RADICAL FORMED BY HYDROGEN ATTACHMENT IN GAS PHASE

Daiki Asakawa⁽¹⁾ - Hidenori Takahashi⁽²⁾ - Shinichi Iwamoto⁽²⁾ - Koichi Tanaka⁽²⁾ - Katsumasa Kamiya⁽³⁾

⁽¹⁾ AIST, Research Institute for Measurement and Analytical Instrumentation, Tsukuba - ⁽²⁾ Shimadzu Corporation, Koichi Tanaka Mass Spectrometry Research Laboratory, Kyoto - ⁽³⁾ Kanagawa Institute of Technology, Physics, Atsugi

Keywords: Hydrogen atom, N-C α bond cleavage, fragmentation dynamics. Ab initio molecular dynamics

Introduction

The novel fragmentation technique involving the interaction between peptide ion and hydrogen atom in gas phase, named HAD, has been recently developed. HAD is initiated by hydrogen radical attachment to the peptide ion in the ion trap, and the fragment ions due to N-C α bond cleavage were selectively observed. The fragmentation process of HAD is not yet fully understood. In this study, the detail processes of HAD were investigated by *ab initio* calculation.

Methods

HAD experiments were performed using a prototype MALDI QIT-TOF mass spectrometer based on the design of AXIMA Resonance equipped H \cdot source. H \cdot was generated by passing H $_2$ gas through a heated tungsten capillary and introduced into the ion trap. The pressure inside the ion trap chamber was maintained below 5×10^{-2} Pa during the reaction. The fragmentation pathway and dynamics were calculated by Gaussian 16 and Car-Parrinello molecular dynamics (CPMD) programs.

Results

To identify the most probable dissociation mechanism of HAD processes, the peptides, AcLys, AcArg, and AcArg-Ala were used as the models. The protonated model peptide was produced by MALDI with CHCA, and then reacted with H \cdot in ion trap. The HAD-MS 2 spectrum of [M+H] $^+$ showed fragment ions due to N-C α bond cleavage. The N-C α bond cleavage was initiated by H \cdot association to carbonyl oxygen in protonated peptide. The TS barrier for the H \cdot attachment was in the range of 0.45-0.6 eV. Since the temperature of H \cdot introduced to ion trap was estimated to 2300 K, 10% of H \cdot in ion trap have the kinetic energy higher than 0.6 eV, which can react with protonated peptides, leading to aminoketyl radical intermediate. To investigate the detail mechanism of aminoketyl radical fragmentation, we performed the simulation based on CPMD. According to the result, the aminoketyl radical immediately underwent homolytic N-C α bond cleavage. The homolytic cleavage resulted in the donation of an electron from the N-C α bond to the nitrogen atom, producing an c' fragment and z \cdot radical.

Conclusions

The use of hot H^{\bullet} efficiently react with protonated peptide in gas phase, inducing N-C α bond cleavage. The dissociation rate on N-C α bond was strongly influenced the temperature of aminoketyl radical. When the temperature was set to be > 500 K, the N-C α bond occurred within 1 ps. Because the heading of ion trap enhanced yield of fragmentation, the computational results were in good agreement with the HAD-MS² result of protonated peptides.

Novel Aspect

H^{\bullet} attachment to peptide leads to aminoketyl radical, which immediately underwent homolytic bond cleavage at N-C α bond

WOr.08-3 ONE-ELECTRON REDUCTION OF ISOLATED METALLO-SUPRAMOLECULAR AGGREGATES

Marianne Engeser

Universität Bonn, Kekulé-Institut für Organische Chemie und Biochemie, Bonn

Keywords: metallo-supramolecular aggregates, supramolecular chemistry, ESI, ECD, ETD

Introduction

Are the metal centres in metallo-supramolecular aggregates electronically isolated from each other? Is there an information transfer via ligands or space? How redox-innocent are the ligands? Gaining knowledge on the redox behaviour of metallo-supramolecular aggregates in solution is highly desirable, but rather difficult due to the inherent lability of the system: It is hard to oxidise or reduce only the aggregate of interest – and not affect one of the other compounds of the (desired) manifold of species present in equilibrium. Herein, a means is presented to study the reduction of isolated supramolecular aggregates under well-defined conditions using mass-spectrometry.

Methods

High-resolution electrospray ionisation (ESI) mass spectra were recorded on a Bruker Apex IV FT-ICR or a Thermo Fisher Scientific LTQ Orbitrap XL mass spectrometer. A series of multiply charged cations typically is observed for metallo-supramolecular aggregates due to abstraction of several anions during the ESI process. After mass-selection of the species of interest, one-electron reduction is performed in the gas phase using the electron capture dissociation (ECD) technique in the FT-ICR or the electron transfer dissociation (ETD) method in the ion trap module of the Orbitrap XL instrument. Reduction products of interest were again mass-selected and fragmented by collision-induced dissociation (CID) or subjected to a second one-electron reduction step.

Results

During ETD and ECD, a single electron is attached to the metallo-supramolecular aggregate. The resulting reduced species is either characterised by high resolution tandem mass-spectrometry or can even be subjected to another reducing step (ECD/ECD) prior to analysis. The presented examples include the stepwise reduction of a copper(II) helicate revealing otherwise unaccessible mechanistic details of its interconversion to the respective copper(I) congener with different stoichiometry.[1] Single-electron reduction of homoleptic metallosupramolecular squares shows the non-innocence of bipyridine ligands [2] whereas a similar experiment with bimetallic squares [3] enables the differentiation between different reduction sites. The effect of linker length is illustrated by reduction of dinuclear gold macrocycles of various ring sizes.[4]

Conclusions

Both ECD and ETD can be used to generate highly interesting partly-reduced open-shell metallo-

supramolecular aggregates. The stability and further reactivity of the reduced species in dependence of type of metal building blocks, bridging and blockings and aggregate size and charge are addressed.

Novel Aspect

ECD and ETD typically are applied to characterize peptide samples whereas they are used to generate and further study the chemistry of highly unusual reduced metallo-supramolecular aggregates herein.

References

1. Gütz C., Hovorka R., Struch N., Bunzen J., Meyer-Eppler G., Qu Z.-W., Grimme S., Cetina M., Topiæ F., Rissanen K., Engeser M., Lützen A., *J. Am. Chem. Soc.*, 136, 1830-11838 (2014).
2. Hovorka R., Engeser M., Lützen A., *Int. J. Mass Spectrom.*, 354, 152-158 (2013).
3. Angurell I., Ferrer M., Gutiérrez A., Martínez M., Rocamora M., Rodríguez L., Rossell O., Lorenz Y., Engeser M., *Chem. Eur. J.*, 20, 14473-14487 (2014).
4. Kobialka S., Müller-Tautges M., Schmidt M.T.S., Schnakenburg G., Holloczki O., Kirchner B., Engeser M., *Inorg. Chem.*, 54, 6100-6111 (2015).

WOr.08-4 REACTION DYNAMICS OF IONIZED 1-METHYLPYRENE: WHY IS THE 1-METHYLENEPYRENE ION FORMED FROM H ATOM LOSS AND NOT A TROPYLIUM-CONTAINING ION?

Paul Mayer⁽¹⁾ - **Brandi West**⁽¹⁾ - **Bethany Lowe**

University of Ottawa, University of Ottawa/Dept. of Chemistry and Biomolecular Sciences, Ottawa

Keywords: polycyclic aromatic hydrocarbon ions, astrochemistry, iPEPICO, collision-induced dissociation

Introduction

The chemical origins of the Diffuse Interstellar Bands have remained elusive. Polycyclic aromatic hydrocarbons, PAHs have been postulated to be sources of DIBs due to the similarity between the IR emissions of the cations of PAHs and the DIBs. One molecule that has been investigated for its potential role as a DIB source is 1-methylpyrene. Previous work has suggested that, unlike toluene ions, ionized methylpyrene dissociate by H loss from the methyl group [1-5]. This study aims to find out why.

Methods

iPEPICO experiments were conducted on the VUV beamline at the Swiss Light Source (Paul Scherrer Institut, Villigen, Switzerland). The mass-selected TPES was acquired to establish the IE. RRKM theory was used to determine the 0 K activation energies (E_0) and entropies of activation ($\Delta^\ddagger S_{1,000K}$) by modeling the experimental breakdown diagrams. Reaction pathways were calculated at the UCCSD/6-31G(d)//B3-LYP/6-31G(d) level of theory.

Results

The energetics derived from modelling ($E_0 = 2.84 \pm 0.36$ eV) is the first experimental value reported for the dehydrogenation of ionized 1-methylpyrene. The loss of H could form either 1-methylenepyrene (A) or a tropylium-motif containing ion (B). The transition state of highest energy in the pathway to B lies ~ 3.20 eV, while that for A is only 2.54 eV, favoring simple H loss from the methyl group, in contrast to that observed for ionized toluene (2.12 eV compared to 2.17 eV [6]). The reason for this flip in relative energies for two the precursors is evident when the geometries for the key TS are compared. In the case of ionized toluene, the TS has equivalent bond angles consistent with a planar 7-membered ring (128.7°). The presence of the three additional rings in 1-methylpyrene prevents the ring from expanding, with the most

constrained angle being $\angle 123$ with a value of 125.0° . The presence of geometric constraints prevents the ion from accessing routes of molecular relaxation which would lower the transition state energy.

Conclusions

The derived energetics, combined with computationally determined reaction coordinates, allowed for the determination that the dehydrogenation product of 1-methylpyrene cations is the 1-methylenepyrene cation. The tropylium-motif containing fragment ion requires significantly more energy to be formed due to steric hindrance in key transition states caused by the constraints imposed by the three benzene rings. The $\Delta_f H^\circ$ (1-methylenepyrene cation) was derived to be $945 \pm 31 \text{ kJ mol}^{-1}$.

Novel Aspect

The preferential formation of 1-methylenepyrene fragment ions is explained by calculations of the key transition states in the reaction pathway.

References

1. Leger A, D'Hendecourt L, Defourneau D, *A&A*, 293, L53 (1995).
2. Parisel O, Ellinger Y, *The Diffuse Interstellar Bands: Contributed Papers*. Tielens, A.G.G.M. ed. NASA Conference Publication 10144. 91 (1994).
3. Kokkin DL, Simon A, Marshall C, Bonnamy A, Joblin C, *Proc. Int. Astronomical Union*, 9, 286 (2014).
4. Rapacioli M, Simon A, Marshall CCM, Cuny J, Kokkin D, Spiegelman F, Joblin C, *J. Phys. Chem. A*, 119, 12845 (2015).
5. Jusko P, Simon A, Wenzel G, Brünken S, Schlemmer S, Joblin C, *ArXiv e-prints* [Online] (2018).
6. Lifshitz, C., *Acc. Chem. Res.* 1994, 27, 138

Wor.08-5 HIGHLY CHARGED PROTEIN IONS: THE STRONGEST ACIDS TO DATE

Muhammad Zenaidee ⁽¹⁾ - William Donald

UNSW, School of Chemistry, Sydney

Keywords: Supercharging, ion-molecule reactions, protein sequencing, Electrospray ionisation, Gas-phase reactions

Introduction

In electrospray ionisation (ESI), highly charged protein ions (HCPI) can be formed that fragment readily to result in higher sequence coverage than ions with less charge. Recently, we discovered that HCPIs can be formed in such high charge states [1] that they protonate N_2 , O_2 and Ar at room temperature. [2] Here, we report a low pressure ESI source that can be used to form HCPIs in even higher charge states that protonate He; i.e. HCPIs are the most acidic entities known.

Methods

Ion molecule reactions of HCPI's with atmospheric gases were conducted on a linear quadrupole ion trap mass spectrometer (LTQ-MS; Thermo Scientific) equipped with an external ESI source. A low-pressure electrospray ionisation source (LP-ESI) was developed using computer aided design and machined using a CNC-mill. Calculated gas-phase values were determined by software written in Python and can be obtained at <https://github.com/mgleeming/PredictPrPlus>.

Results

The extent of protein ion charging in ESI is limited by proton-transfer reactions with atmospheric gases (e.g., N_2 and O_2) surrounding the ESI source. [2] ESI at low pressures should improve

protein ion charging by reducing the number of proton-transfer reactions that can occur between protein ions and N₂ and O₂. Based on the data, it suggests that HCPIs are formed via the ion evaporation model, and not the charge residue model as HCPIs are more than 6 eV less basic than the least basic solution additives.

Using the LP-ESI source, protein ion charge states can be increased by over 40%. For carbonic anhydrase, the charge state improved from 44+ to 58+. Highly charged proteins ions formed by low pressure ESI are sufficiently acidic that He can be protonated in room-temperature ion-molecules reactions; i.e., these ions are the most acidic entities that have been isolated to date.

Conclusions

Using the LP-ESI for the formation of extremely high charge states should allow for improvements in top-down mass spectrometry because protein ions in high charge states are highly reactive and fragment readily. Thus, it is anticipated that the formation of HCPIs by low pressure ESI and chemical supercharging will significantly improve top-down protein sequence analysis by tandem mass spectrometry.

Novel Aspect

The limit to protein charging was found to be due to ion molecule reactions between HCPIs and atmospheric gases. Using LP-ESI, protein ions can be formed in higher charge states.

References

1. Zenaidee M., Donald W., *Analyst*, 140, 1894-1905 (2015)
2. Zenaidee M., Leeming M., Zhang F., Funston T., Donald W., *Angew. Chem. Int. Ed.*, 56, 8522-8526 (2017)

WOr.09 - FTMS and innovative methods for (bio)analyses

Chairs: Ljiljana Paša-Tolić, Giuliana Bianco

WOr.09-1 **Keynote: NATIVE TOP-DOWN FOURIER TRANSFORM MASS SPECTROMETRY OF PROTEINS AND PROTEIN COMPLEXES AS A TOOL FOR STRUCTURAL BIOLOGY AND DRUG DEVELOPMENT**

Joseph A. Loo - Huilin Li - Wonhyeuk Jung - Carter Lantz - Michael Nshanian

University of California, Los Angeles (UCLA), Los Angeles, CA USA

Keywords: native mass spectrometry, FT-ICR, protein complexes, ligand binding, membrane proteins

Introduction

Native mass spectrometry (MS) of protein assemblies reveals size and stoichiometry. But elucidating their structures to understand their function is more challenging. Membrane proteins are involved in many biological functions, but their characterization by a variety of techniques is challenging. We use top-down Fourier transform ion cyclotron resonance (FT-ICR) MS to derive ligand-binding sites and topological information for large proteins and complexes.

Methods

Proteins were analyzed by nanoESI-MS under native buffer conditions (pH 6.8, 200 mM ammonium acetate) using a high resolution 15-Tesla Bruker Solarix FT-ICR system. The FT-ICR MS instrument is equipped with electron capture dissociation (ECD; 1-3 eV) and a CO₂-laser operating at 10.6 μm for activated ion-ECD and infrared multiphoton dissociation (IRMPD). Electron ionization dissociation (EID) was effected by using higher energy electrons (30-50 eV).

Results

Native top-down MS generates information on the surface topology, ligand binding sites, and post-translational modifications of protein complexes [1]. We are using native MS/MS to investigate the molecular action of compounds that prevent amyloid fibril formation in neurodegenerative diseases such as Alzheimer's and Parkinson's disease [2]. To date, membrane proteins have not been addressed by native top-down MS. Optimal activation of membrane proteins in detergent micelles is required. For the homotetrameric 97 kDa water channel, Aquaporin-Z (AqpZ), increasing the collision energy yielded varying degrees of 4M liberation from the detergent micelles as well as fragmentation into the asymmetric dissociation 3M and M substituents. Activated ion ECD to release weakly-bound products generated a top-down spectrum rich in c-/z-products. Product ions originate from regions in the C-terminal portion of the structure, with most localized to the periplasmic domains.

Conclusions

Electrospray ionization's gift for transforming solution-phase macromolecules into *gas-phase ionized* counterparts without disrupting weak noncovalent interactions is key for applying MS to study protein complexes. For the first time, native top-down MS yields structural and sequence specificity for characterization of membrane proteins. The various activation methods available with FT-ICR MS allow protein complexes to be measured with high accuracy.

Novel Aspect

Native top-down mass spectrometry yields higher order structural parameters and ligand binding sites for protein-ligand complexes and membrane proteins.

References

1. Li H., Nguyen H.H., Ogorzalek Loo R.R., Campuzano I.D.G., Loo J.A., *Nature Chemistry*, 10, 139-148 (2018).
2. Acharya S., Safaie B.M., Wongkongkathep P., Ivanova M.I., Attar A., Klärner F.-G., Schrader T., Loo J.A., Bitan G., Lapidus L.J., *Journal of Biological Chemistry*, 289, 10727-10737 (2014).

WOr.09-2 ADVANCES IN ORBITRAP INSTRUMENTATION FOR ANALYSIS OF HIGH-MASS PROTEINS AND PROTEIN COMPLEXES

Alexander Makarov⁽¹⁾ - Eugen Damoc⁽¹⁾ - Eduard Denisov⁽¹⁾ - Kyle Fort⁽¹⁾ - Michiel van de Waterbeemd⁽²⁾ - Dmitriy Boll⁽¹⁾ - Dmitry Grinfeld⁽¹⁾ - Maria Reinhardt-Szyba⁽¹⁾ - Mikhail Belov⁽¹⁾ - Albert J. R. Heck⁽²⁾

⁽¹⁾ Thermo Fisher Scientific, Hanna-Kunath Str. 11, Bremen - ⁽²⁾ Biomolecular MS and proteomics, Utrecht Univ, Utrecht

Keywords: Orbitrap analyzer, high resolution, native mass spectrometry, protein complex, top-down

Introduction

Analysis of high-mass protein and protein complexes poses unique challenges to mass spectrometry (MS), driving deep re-thinking of principles earlier validated on small molecules and peptides. In this work we examine limitations of standard MS for high-mass analysis and explore new ways for extending high-mass performance in various directions.

Methods

Experiments were performed using a modified Thermo Scientific™ Q Exactive™ Plus Orbitrap™ MS with a series of enhancements, including pulsed trapping of ions in the front end of the MS, reduction of RF frequencies on all ion guides and changes in Orbitrap injection. Separately, a new geometry of the Orbitrap analyzer with better vacuum within the trap was explored.

Results

Advances towards analysis of ions in the megadalton mass range and m/z detection in excess of 20,000 became possible only after tackling several major technical challenges, especially insufficient or poorly controllable desolvation and low transmission of high-mass ions. These challenges were addressed by instrument modifications listed above and guided by ion-optical modeling.

Performance was evaluated in MS/MS and pseudo-MS3 experiments for several homomeric and heteromeric megadalton non-covalent protein complexes, including ribonucleoprotein ones.

Separately, a new research-only Orbitrap design facilitated better mechanical balance and vacuum within the trap, enabling for the first time baseline isotopic resolution on intact 150 kDa antibodies both under denaturing and native conditions in infusion mode while employing conventional nitrogen supply for C-trap and HCD cell.

Conclusions

Significant improvements in sensitivity and resolving power were demonstrated for intact proteins, MS/MS and pseudo-MS3 for native top-down analysis of membrane proteins and megadalton non-covalent protein and ribonucleoprotein complexes. This progress yields new insights into the composition, dynamics and interactions of some of nature's most elaborate and enigmatic molecular machines.

Novel Aspect

Theoretical and practical aspects of improving Orbitrap performance for high mass analysis.

WOr.09-3 ION TRAPS WITH QUADRATIC POTENTIALS IN FOURIER TRANSFORM BASED MASS SPECTROMETRY. HARMONIZED MULTI-ELECTRODE KINGDON TRAPS AND FT ICR CELL.

Evgeny Nikolaev^(1,2,4) - ***Gleb Vladimirov***⁽¹⁾ - ***Oleg Kharybin***⁽¹⁾ - ***Luis Fernando Velasquez-Garcia***⁽³⁾ - ***Anastasia Fursova***⁽¹⁾ - ***Zhumei Sun***⁽³⁾ - ***Mikhail Sudakov***⁽¹⁾ - ***Anton Lioznov***⁽³⁾

⁽¹⁾ Skolkovo Institute of Science and Technology, Moscow, Russia - ⁽²⁾ Institute for Energy Problems of Chemical Physics, Moscow, Russia - ⁽³⁾ Massachusetts Institute of Technology, Massachusetts Institute of Technology, Cambridge MA - ⁽⁴⁾ Moscow Institute of Physics and Technology, Moscow, Russia

Keywords: FT MS, ion traps, Harmonized Kingdon traps, dynamically harmonized FT ICR cell

Introduction

The review of current stage of development of the ion trap devices, which are using electric fields with quadratic dependence on coordinate and FT MS approach to signal detection and processing, will be given with emphasize on multiple electrode harmonized Kingdon traps[1] and dynamically harmonized FT ICR cells [2]. The general approach both analytical and digital to electric field calculation for such devices will be presented as well.

Methods

Trap electrode geometries were obtained as the result of solving Laplace and Hamilton-Jacobi equations. Regular high precision machine work and three different 3D printing Methods Selective Laser Sintering, Lost-wax Micromolding (Microcasting) and Inkjet Metal Printing as a methods of traps electrode making are discussed. Simulation of electric field distribution and ion motion dynamics was done by using combined SIMION and SIMAX software.

Results

Ion motion in the electric field with quadratic potential is oscillation, which frequency depends only on m/z and the field curvature. By measuring this frequency and applying FT MS processing method, we can determine m/z of oscillating ions. For static electric field it is possible to create quadraticity in only one direction what follows from basic laws of electrodynamics. One

dimensional quadratic electric fields have been offered for the first time by Knight [3]. Golikov [4] and Koester [5] have offered the options for harmonized traps with multiple internal electrodes. The general approach both analytical and digital to electric field calculation for such devices will be presented. Many variants of the four electrode harmonized Kingdon trap were designed and manufactured. Initial evaluations of different traps performances will be given. The idea of electric field quadraticity emerged in FT ICR mass spectrometry through dynamically harmonized FT ICR cell as a method of significant improvement of resolving power will be discussed as well.

Conclusions

Harmonized multi-electrode ion traps of Kingdon type with quadratic dependence of electric potential along one of the coordinates could be manufactured by different modern technologies both conventional and additive with electrode surface precision high enough to deliver resolving power comparable to Orbitrap [6] analyzer. Dynamically harmonized FT ICR cell and harmonized Kingdon traps could be described analytically using common mathematical approach.

Novel Aspect

The new approaches to characterization of FT MS ion traps based on harmonized multi-electrode Kingdon traps and dynamically harmonized FT ICR traps was developed.

References

1. Kingdon K., *Physical Review*. 21 (4), 408–418 (1923).
2. Nikolaev E., Boldin I., Jertz R., Baykut G. *J Am Soc Mass Spectrom*, 22(7), 1125-1129, (2011).
3. Knight, R., *Applied Physics Letters*.38 (4): 221–223, (1981).
4. Nikitina D. (Golikov) *Ion traps in dynamic mass spectrometry Ph.D. Thesis*, Saint Petersburg State Polytechnical University, 165p, (2006).
5. Köster C., *Int J Mass Spectrometry*, 287, 114-118, (2009).
6. Makarov, A. *Analytical Chemistry*. 72 (6): 1156–62, (2000).

WOr.09-4 THE NEXT DIMENSION IN PROTEOMICS

Pui Yiu Yuko Lam - Christopher A. Wootton - Federico Floris - Cookson K. C. Chiu - Tomos Morgan - Rémy Gavard - Maria van Agthoven - Meng Li - Peter B. O'Connor

University of Warwick, Department of Chemistry, Coventry

2-dimensional mass spectrometry (2DMS) offers intriguing new advantages in proteomics analysis. In particular, 2DMS provides the ability to fragment all proteins/peptides simultaneously in a complex mixture and to accurately correlate which fragment derives from which precursor. The resulting mass spectra are somewhat unusual in appearance, resembling the 2D-NMR results from which the experiment is derived, but the spectra are relatively easy to understand in the context of modern MS/MS experiments. In this presentation, we will discuss the current results on how to best optimise 2DMS experiments for proteomics (in both Bottom-up and Top-down modes). In particular, we will explore the comparison between 2DMS-based proteomics and more traditional nLC-MS/MS results.

For this comparison, a series of standards including BSA and mixtures of standard proteins were used so that the proteomics comparison would focus on known proteins and peptides. In addition, yeast is compared between the two methods where the proteins are also largely known.

The current results clearly show that, in bottom-up proteomics mode, 2DMS can provide very similar levels of sequence and cleavage coverage on these standards, albeit at lower sensitivity at the moment. Additionally, 2DMS doesn't suffer from the co-eluting peak problem in LC-MS/MS, or rather 2DMS is a method that is uniquely capable of handling many precursor and

fragment species simultaneously by expanding the spectra into a second 'modulation' dimension which allows correlation of precursors with fragments.

For top-down proteomics, we can use MS/MS in the hexapole collision cell, followed by 2DMS in the ICR cell to effectively achieve MS³ data for all primary fragments for each and every precursor protein. Our current results show that this methodology can improve cleavage coverage for proteins in top-down proteomics experiments, but further research is ongoing. Latest results will be presented.

WO.09-5 **ULTRA-HIGH RESOLUTION ORBITRAP FTMS: IMPLEMENTATION AND APPLICATIONS**

Yury O. Tsybin⁽¹⁾ - **Konstantin O. Nagornov**⁽¹⁾ - **David P. A. Kilgour**⁽²⁾ - **Natalia Gasilova**⁽³⁾ - **Laure Menin**⁽⁴⁾ - **Anton N. Kozhinov**⁽¹⁾

⁽¹⁾ Spectroswiss Sàrl, EPFL Innovation Park, 1015 Lausanne, Switzerland - ⁽²⁾ Nottingham Trent University, School of Science and Technology, Nottingham - ⁽³⁾ Ecole Polytechnique Federale de Lausanne, Mass Spectrometry Service Facility, Sion - ⁽⁴⁾ Ecole Polytechnique Federale de Lausanne, Mass Spectrometry Service Facility, Lausanne

Keywords: FTMS, Orbitrap, ultra-high resolution, top-down, lipidomics

Introduction

Previously, we implemented a high-performance data acquisition technology enabling users to access unprocessed data (transients) from any FTMS instrument (ICR or Orbitrap). Now, we extend this capability to enable acquisition of much longer, than typical, transients on their FTMS instruments. This instrumentation development, together with advanced data processing software, unlocks previously unseen ultra-high resolution (UHR) FTMS performance.

Methods

Drugs, metabolites, complex mixtures, proteins, and tissues were analyzed using commercial ICR and Orbitrap FTMS instruments; including Q Exactive and Fusion Orbitrap FTMS series. Mass spectral and transient data were acquired in parallel using original manufacturer and in-house developed data acquisition hardware (FTMS Booster). Custom software was developed for acquisition and processing of very large datasets of transients and mass spectra.

Results

Orbitrap FTMS applications are normally performed with transients (much) shorter than 1 s. Our preliminary results demonstrate that modern Orbitrap mass analyzer scan sustain the high coherence level of periodic ion motion - enabling routine acquisition of 3-10 s transients, even from benchtop Orbitraps. As a result, resolution levels exceeding 500000 at m/z 800 can be achieved. This ultra-high resolution (UHR) enables isotopic fine structure analysis, up to at least 1000 Da, providing increased molecular identification confidence for, e.g., drugs and metabolites (lipids). This UHR performance may also support quantitative analyses of lipids and peptides subjected to neutron-encoded isotopic labeling that yields nearly isobaric species. Interestingly, top-down protein analysis and imaging at these UHR levels is an uncharted territory, which is now ready for exploration. Among other examples, we will demonstrate increased sequence coverage and confidence in product ion assignment for top-down UHR FTMS using HCD MS/MS on a Q Exactive HF.

Conclusions

High manufacturing quality of Orbitrap mass analyzers together with high performance of modern data acquisition and processing technologies enable routine UHR FTMS for applications beyond proteomics. We will show how the UHR FTMS may leverage life sciences applications, including metabolomics, top-down, and imaging. Interestingly, even the benchtop Orbitrap

technology (for selected mass spectrometers) is capable of delivering the UHR FTMS performance.

Novel Aspect

Implementation, performance evaluation, and application of ultra-high resolution capability on Orbitrap FTMS instruments, including on a benchtop series

WOr.10 - Food additives and contaminants

Chairs: Jose Diana Di Mavungu, Roberta Galarini

WOr.10-1 **Keynote:** DIRECT FOOD ANALYSIS BY (TRANSPORTABLE) AMBIENT IONIZATION MASS SPECTROMETRY

Michel W.F. Nielen⁽¹⁾ - **Marco Blokland**⁽²⁾ - **Arjen Gerssen**⁽²⁾ - **Michael van Dam**⁽²⁾ - **Germán Augusto Gómez-Ríos**⁽³⁾ - **Janusz Pawliszyn**⁽³⁾

⁽¹⁾ Wageningen University, Wageningen, The Netherlands - RIKILT, Wageningen, The Netherlands - ⁽²⁾ RIKILT, Wageningen, The Netherlands - ⁽³⁾ University of Waterloo, Waterloo, Canada

Keywords: food analysis, ambient ionization, mass spectrometry, sample preparation, on-site detection

Introduction

Many health and food related apps have been developed in recent years, followed by the development of wearable sensors for blood pressure, heart beat, moisture and biochemical parameters such as blood glucose levels. Obviously, mankind gets used to measuring all kinds of parameters themselves and subsequently take control over their own life, based on analytical data. In the near future, also the pre-screening of residues and contaminants in food will be performed on-site at the farm, in industry, in retail and, eventually, even done by consumers at home. Typically, three different approaches to selectivity in portable food analytics (each at different development stages), can be distinguished: (1) generic spectral profiles from smartphone-linked optical spectroscopy sensors, such as near infrared (NIR) combined with chemometric models; (2) smartphone-based biorecognition assays; and (3) higher resolution spectral data, such as available from mass spectrometry (MS). With respect to the MS option, there is an urgent need for (trans)portable MS, elimination of tedious sample clean-up, simplified sample introduction, and data evaluation tools that communicate wireless to smartphones, tablets and servers of relevant stakeholders, including the consumer. Sample introduction can, at least in theory, partly be replaced and simplified by one of the >70 Ambient Ionization MS (AIMS) techniques.

Methods

In this research, different AIMS techniques with potential for coupling to (trans)portable mass analyzers have been critically investigated for their applicability towards the analysis of pesticides, banned veterinary drugs such as β -agonists and steroids, regulated drugs such as analgesics and antibiotics, and natural toxins in various food matrices. Four different AIMS techniques were explored: Matrix Assisted Ionization (MAI), Handheld Desorption Atmospheric Pressure Chemical Ionization (DAPCI), ID-CUBE Direct Analysis in Real Time (DART) and Coated Blade Spray (CBS), either in combination with a benchtop Orbitrap or a transportable single quadrupole MS.

Results

MAI is one of the most simple ionization options since no high voltage, lasers or gas supplies are

needed, the sample or extract is mixed with a matrix and simply tapped against the opening of the cone or heated inlet of the MS. Secondly, a handheld battery powered DAPCI device, just using ambient air, can be used to ionize analytes directly from a sample surface. Third, a mini-DART setup (a modified version of the ID-CUBE) can be used but would require a helium gas supply. With this technique the sample is simply transferred onto a metal mesh sample card having the size of a credit card. Fourth, the recently developed CBS technique, only requires a high voltage. Both the metal mesh card and the CBS technique offer the benefit of trace enrichment following coating with a selective adsorbent. In this presentation an overview is given about the requirements, challenges and preliminary results of direct food analysis by these different (transportable) AIMS approaches.

Conclusions

All tested devices were able to successfully ionize some of the tested analytes, some of them even in various food commodities. A general drawback of these ambient ionization techniques is the poor reproducibility, but, where feasible, isotope dilution of the sample may compensate for that to a large extent. Detection limits will be typically (much) higher when compared to state-of-the-art high-end MS instruments so applications may be limited to regulated substances having relatively high maximum limits values, unless selective enrichment approaches, such as CBS, are being exploited. Nevertheless, the initial applications underline the great potential of ambient ionization with (trans)portable MS. Hence, it is realistic to envisage that portable MS, already successful in battlefield- and forensic applications for many years, will soon show up in food control as well.

Novel Aspect

A paradigm shift in food control, from the laboratory to on-site pre-screening using, among others, (trans)portable MS is upcoming.

WOr.10-2 GUM ARABIC CHARACTERISATION AND ITS DETECTION IN WINE USING HIGH RESOLUTION MASS SPECTROMETRY

Tiziana Nardin - Giulia Ferrari - Roberto Larcher
Fondazione E. Mach, CTT, San Michele a/A

Keywords: adjuvant, Orbitrap™, characterization, wine addition

Introduction

Gum Arabic is a dried exudate obtained from several *Acacia* species from the sub-Saharan region of Africa [1]. It is a complex mixture of macromolecules, mainly carbohydrates and proteins, widely used as an edible ingredient in various types of food as well as for non-food applications [2]. For oenology, OIV approves the use of both *A. senegal* and *A. seyal* gums. Europe indicates a maximum of 200 mg/L as the recommended technological dosage [3].

Methods

The study was performed using a Thermo Ultimate R3000 UHPLC equipped with a byphenyl column (3 x 150 mm, 2.7 µm). Separation was obtained at a flow rate of 0.3 mL/min with a ternary mobile phase with 2% formic acid, acetonitrile and water. Mass spectra were acquired through a full MS experiment at 70,000 FWHM resolution using a Q-Exactive™ HRMS equipped with a heated electrospray ionization (HESI-II) interface working in both positive and negative ion modes.

Results

The full scan MS profile, both in positive and negative ionization, of 45 gum Arabic samples was evaluated in order to check if specific ions could selectively characterise the botanical (*A. seyal*

and *A. senegal*) and/or the geographic origins of gums (Kordofan, Hashab and Vereke). The PLS-DA statistical approach indicated 18 masses as possible predictive origin markers, allowing the correct reclassification of 98% of the samples to the corresponding cases. Moreover, we analysed 40 micro-winery genuine wines (8 varieties) and then reanalysed the same ones after randomly adding 200 mg/L of gum Arabic of the 2 origins (*seyal*, n=10; *senegal*, 10). Comparing their ion profiles, 5 of the previous 18 masses, not detected in genuine wines, permitted the identification of the addition of gum Arabic. PCA based on the 5 masses allowed us to distinguish between wines added with *A. senegal* and *A. seyal*. Finally, each isotopic patterns and ion fragmentation allowed us to define their molecular formulas and structures.

Conclusions

The HRMS untargeted approach allowed us to identify, through PLS-DA analysis, 18 masses (*m/z*) that could be proposed as possible markers for the differentiation of the commercial gum Arabic samples, both with respect to the botanical (*A. seyal* and *A. senegal*) and geographic (Kordofan, Hashab and Vereke) origin. 5 of these 18 masses (*m/z* 152.1067, 166.1225, 181.0459, 585.1428, 643.1403) were also effective in detecting the presence of 200 mg/L gum Arabic in wine.

Novel Aspect

A new rapid untargeted high resolution MS approach allows us to predict gum Arabic origin and detect its presence in wine.

References

1. Sprenger S., Meylahn K., Zaar A., Dietrich H., Will F., *Eur Food Res Technol*, 240, 909-921 (2015).
2. Coppen J.J.W., *Food and Agriculture Organization of the United Nations*, Rome, Italy (1995).
3. *OIV-3.Codex*, Parigi, Francia (2018).

WOr.10-3 THE FRAGMENTATION BEHAVIOUR OF PYRROLIZIDINE ALKALOIDS AND ITS IMPLICATIONS ON FOOD ANALYSIS

Christoph Czerwenka

Austrian Agency for Health and Food Safety, Group for Contaminant and Special Analysis, Wien

Keywords: pyrrolizidine alkaloids, fragmentation, isomers, food analysis

Introduction

Pyrrolizidine alkaloids are natural toxins which are produced by numerous plant species and may contaminate food such as honey or herbal teas. This substance group exhibits a huge structural variety, which, however, can be clustered into several classes which share common structural elements. Several of these alkaloids are isomeric to one another, thus complicating their differentiation during food analysis.

Methods

The analysis of pyrrolizidine alkaloids is commonly performed by LC-MS/MS. Therefore, the in-depth knowledge of their fragmentation behaviour is of great importance for a selective and sensitive analysis. The fragmentation patterns of more than 30 pyrrolizidine alkaloids were studied by low- and high-resolution mass spectrometry on a triple-quadrupole and a QTOF instrument. Food samples were analysed by LC-MS/MS utilising the previously obtained information.

Results

The obtained product ion spectra and accurate masses of the fragments of the different pyrrolizidine alkaloid species allowed the deduction of the structures of the fragments and their

assignment as characteristic product ions of certain pyrrolizidine alkaloid classes. These structural relationships will be presented and discussed in detail. It was observed that isomeric pyrrolizidine alkaloids – both stereoisomers as well as structural isomers – have all the same fragments; however, the relative abundancies of the different fragments were found to vary dramatically between the isomers. For a correct differentiation between these isomers in the course of the LC-MS/MS analysis a chromatographic baseline separation is essential, otherwise the obtained quantitative results will strongly depend on the ratio of the isomers and the analytes may fail common identification criteria based on MRM ratios. Results of real-life honey and herbal tea samples will be shown that highlight the need for a careful differentiation between the isomers.

Conclusions

A detailed understanding of the fragmentation of pyrrolizidine alkaloids is of great value for a correct analysis of food products for these compounds. Stereoisomeric pyrrolizidine alkaloid species can exhibit significantly different mass spectrometric behaviours, thus calling for a sophisticated chromatographic separation to allow a careful differentiation between them in order to obtain a true picture of contamination levels in food products.

Novel Aspect

The fragmentation behaviour of pyrrolizidine alkaloids was elucidated in detail and its impact on the correct determination of contamination levels in food demonstrated.

WOr.10-5 THE SWEET DEAL OF ADDING ION MOBILITY TO FULL SCAN ACQUISITIONS TO PROFILE STEVIOL GLYCOSIDES KNOWN AND KNOWN-UNKNOWN IN FOOD MATRICES

S  verine Gosciny⁽¹⁾ - M  lanie Herman⁽¹⁾ - Els Van Hoeck⁽¹⁾ - Mike McCullagh⁽²⁾

⁽¹⁾ Sciensano, Organic Contaminants and Additives, Brussels - ⁽²⁾ Waters Corporation, Science Operations, Wilmslow

Keywords: Steviol glycosides, isomers, Collision Cross Section (CCS), ion mobility spectrometry (IMS), known-unknowns

Introduction

In 2011, the EU authorised steviol glycosides as a sweetener food additive (E960), its natural origin is key to success in the food market. Like many botanical extracts, E960 is a complex mixture of compounds, with diverse physico-chemical/organoleptic properties and is a challenge to analyse. Hence its characterisation in complex food matrices requires a variety of analytical tools. Here, we present a simple screening method using ion mobility mass spectrometry.

Methods

For the assays, RP-C18-Ion Mobility-HRMS analyses were performed on a Synapt G2-Si using standard solutions and real samples. The sample treatment comprises the extraction of the steviol glycosides, followed by C18-SPE clean-up prior to analysis. Mass spectra were collected from m/z 50 to 2000 (resolution 20000 FWHM). The mobility drift gas used was nitrogen and mobility cell calibration was carried out using poly-DL-alanine.

Results

From the analysis of steviol +10 steviol glycosides standards, CCS values, Rt, and fragmentation patterns were derived for precursor, formiate adducts and labile fragments ions. These results were entered into a screening library. Regardless of the molecular complexity of the food additive, co-eluting critical pairs were ion mobility resolved. Using the method, E 960 was profiled in food matrices. The results show the efficiency of CCS values in simplifying the

screening workflow, while other isomeric steviol glycosides were detected. Characterization of these known-unknowns without standards was performed using IMS. Drift time aligned fragment spectra, confirmed the origin of the fragment ions and in combination with Zimmermann [1] fragmentation rules, sugar moiety positions were identified. These compounds were inserted into the library and used for retrospective analysis. The CCS specificity of known-unknowns enabled the distinction of multiple isomers within complex mixtures of food extracts, providing enhanced profiles.

Conclusions

Using IM-HRMS analysis methodology, multiple identification parameters for the same steviol glycoside have been determined. As a result, the screening workflow has been drastically simplified and successfully applied to the analysis of a variety of food matrices containing the E 960 food additive. The additional CCS metric has enhanced confidence in identification of steviol glycosides at trace levels, which is critical considering the sweetness potency of E 960.

Novel Aspect

This work demonstrated the benefits of combining ion mobility to full scan acquisitions to build a comprehensive screening analysis, even when standard references are unavailable.

References

1. Zimmermann B.F., *Rapid Commun. Mass Spectrom.*, 25, 1575–1582 (2011).

THURSDAY, August 30th, 2018

PL05: Plenary Lecture: Thomson Medalist

Chair: Renato Zenobi

PL05 **DRIVING INNOVATION - FROM A PROTEIN SEQUENCE TO A PROTEOME**

John R. Yates

The Scripps Research Institute, La Jolla

Keywords: proteins, proteome, tandem mass spectrometry, informatics

Proteins are the workhorses of the cell. They are the molecules that catalyze the reactions for the metabolic and catabolic activities of the cell, transmit signals and transport necessary molecules and metabolites. Proteins operate in systems, which can be pathways, organelles, cells or even tissues. The development of proteomic technologies allows the study of proteins in the context of cellular systems. This lecture will trace the development of some of the methods of proteomics and demonstrate their use to understand mechanisms underlying the disease Cystic Fibrosis.

ThOr.01 - Glycomics, saccharides and glycoconjugates

Chairs: Catherine Costello, Joseph Banoub

ThOr.01-1 Keynote: HYBRIDATION OF MASS SPECTROMETRY AND ION SPECTROSCOPY FOR GLYCOANALYTICS

Isabelle Compagnon

Institut Lumière Matière. Université Lyon 1/CNRS. Villeurbanne. France. Institut Universitaire de France. Paris

Keywords: Mass Spectrometry; Infrared Ion Spectroscopy; Carbohydrate structure; Carbohydrate sequencing

Deciphering the sequence of carbohydrates is a major challenge in glycosciences today. Whereas spectroscopic approaches (NMR, Raman, ...) and Mass Spectrometry (MS) offer valuable information on carbohydrate structure, they also suffer from severe technical limitations: on one hand, spectroscopy provides refined structural detail but requires a relatively large amount of purified sample, on the other hand MS applies to complex samples in small amounts, but does not disambiguate all carbohydrate isomerisms.

In this context, we have developed a unique hybrid Glycoanalytical approach, which adds a spectroscopic dimension to MS analysis. Based on this technology, it is possible to measure simultaneously the MS spectrum and the spectroscopic fingerprint of a carbohydrate. This enables the disambiguation of all kinds of carbohydrate isomerisms (monosaccharide content, regiochemistry and stereochemistry of the glycosidic bond, and positional isomers of functional modifications) without the need for chemical modification, labeling, or time-consuming purification.

While this new carbohydrate metric mitigates the lack of structural specificity of MS without the limitation of traditional spectroscopy, the identification of a carbohydrate however still relies on the availability of referenced standards. In order to enable true “de novo” analysis, we have established a top-down sequencing strategy based on the demonstration of a memory of the carbohydrate sequence within MS fragmentation products. [1] The performance of de novo carbohydrate sequencing will be illustrated on various classes of carbohydrates.

References

1. B. Schindler et al. Nature Communications, 8 (2017).

ThOr.01-2 APPLICATION OF A NOVEL FUNCTIONALIZED RESIN WITH PHENYLBORONIC ACID FOR CAPTURING OF GLYCATED AND GLYCOSYLATED PEPTIDES AND THEIR MS IDENTIFICATION.

Magdalena Wierzbicka⁽¹⁾ - Monika Kijewska⁽¹⁾ - Mateusz Waliczek⁽¹⁾ - Francesca Nuti⁽²⁾ - Anna Maria Papini⁽²⁾ - Piotr Stefanowicz⁽¹⁾ - Zbigniew Szewczuk⁽¹⁾

⁽¹⁾ Faculty of Chemistry, University of Wrocław, Wrocław - ⁽²⁾ French-Italian Interdepartmental Laboratory of Peptide and Protein Chemistry and Biology, University of Florence, Florence

Keywords: glycopeptides, glycosylation, capturing, PTMs, LC-MS

Introduction

Proteins post-translational modifications, e.g. glycation and glycosylation are investigated as specific biomarkers for diagnosis. Limitations of standard bottom up approach are connected with peptides low relative abundance, low ionization efficiency and sample complexity. The further isolation and concentration methods are required, which can be achieved by phenylboronic acid affinity materials[1].

Methods

The ChemMatrix® Rink Resin functionalization with 4-carboxyphenylboronic acid moiety was performed according to the solid phase peptide synthesis procedure. The glycoconjugates were captured by incubation with the functionalized resin. Capturing and liberation progress were

monitored both by the mass spectrometric Methods LC-MS, FT-ICR-MS, MS/MS (all with ESI ion source) and UV spectroscopy.

Results

Solid supports were efficiently functionalized with various spacers and 4-carboxyphenylboronic acid. The linker was characterized by LC-MS and UV and the substitution level of the resin was determined. The resin is highly stable to various reaction conditions as confirmed by LC-MS results. Synthetic glycoconjugates were synthesized using sugar-modified building blocks [2-4]. The labeling of glycopeptides by dabcyll moiety allowed the monitoring of capturing by UV. The glycosylated peptides were captured efficiently by phenylboronic acid solid support (70%). The capturing progress depends on glycoconjugates concentrations as monitored by UV. The captured products were identified by LC-MS, FT-ICR-MS and MS/MS.

Conclusions

We developed and tested the new method of glycopeptides enrichment. The affinity of glycopeptide to resin functionalized with 4-carboxyphenylboronic acid depends on the structure of saccharide moiety. LC-MS allows identification and quantification of captured of specific glycoconjugates. This method has potential applications in biomarkers identification in medical diagnosis.

Novel Aspect

The application of ChemMatrix® resin allows convenient and fast preparation of supports with various structures, properties and specificity for glycoconjugates capturing.

References

1. Liu Z., He H., *Acc. Chem. Res.*, 50, 2185–2193 (2017)
2. Stefanowicz P., Kijewska M., Kapczyńska K., Szewczuk Z., *Amino Acids.*,38(3):881-889 (2010)
3. Nuti F., Peroni E., Real-Fernandez F., Bonache M. A., Le Chevalier-Isaad A., Chelli M., Lubin-Germain N., Uziel J., Rovero P., Lolli F., Papini A.M., *Biopolymers*, 94 (6), 791-799 (2010)
4. Carganico S, Rovero P, Halperin JA, Papini AM, Chorev M., *J Pept Sci.* 15(2), 67-71 (2009)

ThOr.01-3 GLYCAN ISOMER DIFFERENTIATION USING A SIALIC ACID DERIVATIZATION TOOLBOX FOR HIGHER THROUGHPUT MASS SPECTROMETRIC GLYCOMICS

Noortje de Haan - Karli Reiding - Stefanie Holst - Bram Heijs - David Falck - Guinevere Kammeijer - Viktoria Dotz - Manfred Wuhrer

Leiden University Medical Center, Center for Proteomics and Metabolomics, Leiden

Keywords: sialic acid, ethyl esterification, dimethylamidation, lactonization, isomer differentiation

Introduction

Glycosylation is a highly complex post-translational modification that plays a key role in the organization of multicellular life. Many human proteins are glycosylated, and changes in their glycosylation occur with many disease, a phenomenon which is still poorly understood. We here present mass spectrometric methods for analyzing glycosylation at higher throughput with differentiation of biologically important sialic acid linkage isomers.

Methods

Glycans were enzymatically released from proteins and subjected to sialic acid derivatization in a linkage-specific manner introducing a mass difference for linkage differentiation. The approach was hyphenated with MALDI-MS [1], and -imaging MS [2]. Furthermore, the method was

combined with capillary electrophoresis(CE) – sheathless nanoESI-MS after reducing-end labeling with a positively-charged tag, achieving glycan analysis at ultrahigh sensitivity.

Results

High-throughput glycan analysis after linkage- specific sialic acid derivatization was applied to human serum and plasma samples of thousands of individuals to define disease-associated glycosylation signatures of major human diseases including inflammatory bowel disease, colorectal cancer, rheumatoid arthritis, and diabetes type 2. In addition, glycosylation signatures of age, sex, pregnancy, obesity and inflammation were determined. The derivatization method was transferred to tissue sections in order to allow glycan mass spectrometry imaging with a focus on sialic acid linkage differentiation providing valuable insights into the glycosylation heterogeneity of various types of tumors. The method was further expanded to a CE-MS variant allowed the analysis of the human plasma N-glycome with unprecedented analytical depth which proved particularly valuable for analyzing scarce biological samples that are available in very restricted amounts.

Conclusions

Glycan isomers differing in sialic acid linkages can be readily differentiated using the tailored derivatization toolbox in combination with MS detection. The method revealed glycosylation signatures of major human diseases and provided insights into tumor heterogeneity with respect to glycosylation.

Novel Aspect

A high-specificity derivatization step in the glycan MS workflow allows the rapid differentiation of sialic acid linkage isomers providing novel biomedical insights.

References

1. Reiding KR, Blank D, Kuijper DM, Deelder AM, Wuhrer M. *Anal. Chem.* 86, 5784 (2014)
2. Holst S, Heijs B, de Haan N, van Zeijl RJ, Briaire-de Bruijn IH, van Pelt GW, Mehta AS, Angel PM, Mesker WE, Tollenaar RA, Drake RR, Bovée JV, McDonnell LA, Wuhrer M. *Anal. Chem.* 88, 5904 (2016)

ThOr.01-4 **HILIC MS-HCD/FT-IT AND HYBRID SEARCH ENHANCE THE STRUCTURE ELUCIDATION OF OLIGOSACCHARIDES IN NIST HUMAN MILK REFERENCE MATERIAL**

Connie A. Remorosa - Tytus D. Mak - Stephen E. Stein

Mass Spectrometry Data Center, Biomolecular Measurement Division, National Institute of Standards and Technology, Gaithersburg, MD, 20899-8362, USA

Keywords: Human Milk Oligosaccharides, Fourier Transform-ITMS, HILIC-MS-HCD, Isomers, MS library

Introduction

Elucidation of complex mixture of underivatized isomeric glycans is challenging when commercial standards are not available. The NIST SRM* 1953 consist of human milk oligosaccharides (HMO) from 100 healthy breastfeeding mothers in the US was used. However, the information on these functional molecules in the reference materials is limited. This study focuses on the development of MS database of annotated spectra of oligosaccharides in human milk.

Methods

The human milk SRM 1953 underwent sample preparation [1]. The purified oligosaccharides

were analyzed on Hydrophilic Liquid Interaction Chromatography-electrospray ionization in Fusion Lumos Orbitrap-based MS using HCD&Fourier transform ion trap (FT-IT) fragmentation techniques. The spectra are extracted from LC-MS raw data and processed using in-house algorithm then analyzed using different tools including NIST Tandem mass spectral library and hybrid search.

Results

The high mass accuracy both for precursor and product ions from FT-IT and HCD data was sufficient to enable the assignment of HMOs including novel oligosaccharide structures with different adduct types in negative/positive detection mode. Core group (Gal-Glc) with different terminal sugars, Fuc-Gal-Glc eluted prior to Sia-Gal-Glc based on the HILIC chromatogram. The elution profile of oligosaccharides (>4 sugar units) with 3-4 isomers, such as fucolacto-N-octaoses and sialylated fucolacto-N-hexaoses displayed different retention times on HILIC that enhances oligosaccharide identification. A clear distinction in the HILIC elution pattern of isomeric oligosaccharides related to size and polarity also confirmed the identifications. Moreover, the hybrid search can reliably find modified oligosaccharides in comparison to non-reduced analogs, and vice versa, enabling an increase in match factors.

Conclusions

The HILIC MS/MS and NIST hybrid search facilitates the structural assignment of 70 complex oligosaccharides of 470 spectra including newly identified glycan isomeric structures. The MS library of HMOs demonstrated to aid in identifying glycans in other biological samples, which represents a new variety of mass spectral library – a concept that can be extended to other varieties of biological fluids such as detection of glycans and non digestible artificial sugars in human & bovine milk samples.

Novel Aspect

A comprehensive structural and tandem MS library of annotated spectra of oligosaccharides derived from the NIST human milk reference material: a robust method of detecting other compounds in milk products.

References

1. Wu, S. et al. (2010). Journal of Proteome Research 9, 4138-4151

* SRM (Standard Reference Material) 1953 (human milk); SRM 1849 Infant formula

ThOr.01-5 **COMPETITIVE UNIVERSAL PROXY RECEPTOR ASSAY (CUPRA) FOR QUANTITATIVE HIGH-THROUGHPUT GLYCAN LIBRARY SCREENING**

John Klassen - Elena Kitova - Pavel Kitov - Ling Han - Zhixiong Li

University of Alberta, Chemistry, Edmonton

Keywords: Glycan recognition; library screening; affinity; specificity; ESI-MS

Introduction

The absence of quantitative high-throughput (Q-HT) glycan screening assays is a significant impediment in glycomics. Glycan microarrays, wherein oligosaccharides (OS) are immobilized on glass surfaces, are the dominant HT screening technology. However, glycan microarray screening has many well-documented limitations. Here, we introduce a Q-HT electrospray ionization mass spectrometry (ESI-MS)-based method for OS library screening.

Methods

The *Competitive Universal Proxy Receptor Assay* (CUPRA) assay is based on competitive protein

binding and employs a universal proxy protein ($^{Uni}P_{proxy}$), which binds all components of a library of modified OS (OS^{mod}) through a high affinity tag introduced to the OS. Changes in abundances of the ($^{Uni}P_{proxy}+OS^{mod}$) complexes in solution upon introduction of a glycan-binding protein (GBP) allows for the simultaneous identification and quantification of OS^{mod} ligands.

Results

Implementation of CUPRA was demonstrated using a library of 70 OS^{mod} , produced from OS with structures found in humans, and GBPs with known OS binding specificities. Application of CUPRA to cholera toxin B subunit homopentamer (CTB₅) correctly identified the OS^{mod} containing the Gal β 1-3GalNAc β 1-4(Neu5Ac α 2-3)Gal β 1-4Glc (GM1) and Fuc α 1-2Gal β 1-3GalNAc β 1-4(Neu5Ac α 2-3)Gal β 1-4Glc (Fuc-GM1) structures as the only high affinity ligands in the library. Moreover, the measured affinities agree, within a factor of 2, with the affinities measured for the corresponding oligosaccharides. Application to *Maackia Amurensis* agglutinin (MAA), which is a Neu5Ac α 2-3Gal β 1-4GlcNAc/Glc specific GBP, identified Neu5Ac α 2-3Gal β 1-4Glc (3SL) and Gal β 1-4(Neu5Ac α 2-3)Gal β 1-4Glc (GM2) as ligands. A single ligand, containing the Neu5Ac α 2-6Gal β 1-4Glc structure was found for *Sambucus Nigra* lectin (SNA). This result is consistent with the reported specificity of the GBP for Neu5Ac α 2-6Gal structures.

Conclusions

Together, these results demonstrate that CUPRA allows for the simultaneous screening and quantification of glycan ligands with affinities $\geq 10^3 M^{-1}$. Future efforts will focus on expanding the size and diversity of the OS^{mod} library, with the goal of producing a library of ~ 200 components derived from purified OS. Natural libraries of OS^{mod} will also be produced from OS extracted from human milk and N-glycans released from isolated glycoproteins and tissue. The OS^{mod} libraries will then be used for the discovery of biologically-relevant glycan interactions involving human, bacterial and viral GBPs.

Novel Aspect

Development of the first high-throughput and quantitative glycan screening assay

ThOr.02 - Probing Protein Structure and Dynamics

Chairs: Lars Konermann, Giancarlo la Marca

ThOr.02-1 **Keynote: PROBING THE STRUCTURES OF NATIVE AND PATHOLOGICAL A-SYNUCLEIN SPECIES BY HYDROGEN/DEUTERIUM EXCHANGE MASS SPECTROMETRY**

Thomas J. D. Jørgensen⁽¹⁾ - **Jeppe Buur Madsen**⁽¹⁾ - **Anne Louise Hemdrup**⁽²⁾ - **Cagla Sahin**⁽³⁾ - **Daniel Otzen**⁽³⁾ - **Poul Henning Jensen**⁽⁴⁾

⁽¹⁾ University of Southern Denmark, Department of Biochemistry and Molecular Biology, Odense M, Odense -

⁽²⁾ Aarhus University, Interdisciplinary Nanoscience Center, Aarhus C, Aarhus - ⁽³⁾ Aarhus University, Department of Molecular Biology and Genetics, Aarhus C, Aarhus - ⁽⁴⁾ Aarhus University, DANDRITE - Danish Research Institute of Translational Neuroscience & Department of Biomedicine, Aarhus C, Aarhus

Keywords: intrinsically disordered protein, alpha-synuclein, in-cell, H/D exchange, protein structural dynamics

Introduction

It is well established that purified monomeric α -synuclein is an intrinsically disordered protein, but its structure in the cellular environment has been a matter of intense scientific controversy [1-3] fueled by reports that α -synuclein adopts a stable folded tetramer in cells, but only when

they are intact.[4-7] This controversy is mainly due to the inherent difficulties associated with probing protein structures inside the cell.

Methods

We here present a novel approach to probe the intracellular structure of endogenous proteins and we address the question whether the cellular structure of endogenous α -synuclein is folded or disordered. Our approach is based on in-cell hydrogen/deuterium exchange where α -synuclein is labeled in its cellular environment with D_2O followed by immunoaffinity purification at acidic conditions (in-cell-HDX-IPAC) where the isotopic exchange is quenched.

Results

Recombinant ^{15}N -labeled monomeric α -synuclein was added to human erythrocytes prior to isotopic exchange as an intrinsically disordered extracellular reference protein. This allows for sensitive detection of protection of intracellular endogenous α -synuclein caused by any putative folded stable structure. To ensure that the equilibration of D_2O across the plasma membrane of erythrocytes occurs sufficiently rapid, we determined the deuterium uptake of intracellular and extracellular hemoglobin and it was found to be identical. To enable in-cell HDX-MS analysis of endogenous α -synuclein, an effective enrichment of the protein at quench conditions is required. For this purpose, we developed a novel immunoaffinity approach (in-cell-HDX-IPAC) that maintains strong capture efficiency at acidic conditions (pH 3.5) to limit the amount of back exchange. The enriched α -synuclein was eluted at pH 2.0. Our results show that the deuterium uptake of intracellular and extracellular α -synuclein is identical.

Conclusions

The identical deuterium uptake of intracellular and extracellular α -synuclein suggests that the vast majority of endogenous α -synuclein in erythrocytes is intrinsically disordered. Furthermore, we have also shown that our in-cell-HDX-IPAC approach is able to detect the protection pattern of the presence of nonnative α -synuclein oligomers.[8,9] In addition, we have investigated the stability and dynamics of nonnative α -synuclein oligomers.

Novel Aspect

In-cell H/D exchange followed by immunoaffinity purification at acidic conditions (in-cell-HDX-IPAC) is utilized to probe the intracellular structure of α -synuclein.

References

1. Binolfi A., Theillet F. X., Selenko P. *Biochem. Soc. Trans.*, 40, 950-954,(2012).
2. Theillet F. X., Binolfi A., Bekei B., Martorana A., Rose H. M., Stuijver M., Verzini S., Lorenz D., van Rossum M., Goldfarb D., Selenko P. *Nature*, 530, 45-50 (2016).
3. Fauvet B., Mbefo M. K., Fares M. B., Desobry C., Michael S., Ardah M. T., Tsika E., Coune P., Prudent M., Lion N., Eliezer D., Moore D. J., Schneider B., Aebischer P., El-Agnaf O. M., Masliah E., Lashuel H. A. *J Biol. Chem.*, 287, 15345-15364 (2012).
4. Bartels T., Choi J. G., Selkoe D. J. *Nature*, 477, 107-110 (2011).
5. Dettmer U., Newman A. J., Soldner F., Luth E. S., Kim N. C., von Saucken V. E., Sanderson J. B., Jaenisch R., Bartels T., Selkoe D. *Nat Commun*, 6, 7314 (2015).
6. Dettmer U., Newman A. J., vonSaucken V. E., Bartels T., Selkoe D. *Proc Natl AcadSci U S A*, 112, 9596-9601 (2015).
7. Wang W., Perovic I., Chittuluru J., Kaganovich A., Nguyen L. T., Liao J., Auclair J. R., Johnson D., Landeru A., Simorellis A. K., Ju S., Cookson M. R., Asturias F. J., Agar J. N., Webb B. N., Kang C., Ringe D., Petsko G. A., Pochapsky T. C., Hoang Q. Q. *Proc Natl AcadSci U S A*, 108, 17797-17802 (2011).
8. Mysling S., Betzer C., Jensen P. H., Jørgensen T. J. D. *Biochemistry*, 52, 9097-9103 (2013).
9. Paslawski W., Mysling S., Thomsen K., Jørgensen T.J.D., Otzen D. E. *AngewChemInt Ed Engl*, 53, 7560-7563 (2014).

ThOr.02-2 HDX-MS REVEALS CONFORMATIONAL CHANGES OF IL-23 IN COMPLEX WITH IL-23R/IL-12RB1 AND PHAGE DISPLAY PEPTIDES

Howard Broughton⁽¹⁾ - Leticia Cano⁽¹⁾ - Jeffrey Dodge⁽²⁾ - Isabel Gonzalez⁽²⁾ - Chalmers Michael⁽²⁾ - Feiyu Zhang⁽³⁾ - Sepideh Afshar⁽³⁾ - Alfonso Espada⁽¹⁾

⁽¹⁾ Lilly, S.A., Centro de Investigación, Alcobendas - ⁽²⁾ Eli Lilly and Company, Lilly Research Laboratories, Indianapolis - ⁽³⁾ Eli Lilly and Company, Lilly Biotechnology Center, San Diego

Keywords: Interleukin IL-23, HDX-MS, computational approaches, binding interactions, phage display peptides

Introduction

Interleukin IL-23 is a pro-inflammatory cytokine that contributes to the pathogenesis of autoimmune inflammatory diseases. The inhibition of the activity of IL-23 has promising therapeutic effects for psoriasis, psoriatic arthritis, inflammatory bowel disease, rheumatoid arthritis and multiple sclerosis.[1] IL-23 is a heterodimer composed of the cytokine subunits p19 and p40 that bind to the receptor complex IL-23R and IL-12R β 1 respectively. The discovery of monoclonal antibodies designed to target IL-23 (either p19 and /or p40) have shown significant efficacy in human studies for some of these diseases. Here we present the application of HDX-MS [2] and computational approaches to interrogate the binding interactions of IL-23 with its receptors and the binding epitope of new IL-23 inhibitors such as phage display peptides [3].

Methods

The exchange reaction was performed with a Tecan liquid handler modified to allow for freezing of samples into Millitubes arrayed in a 96 well plate format. Thus, samples were prepared triplicate of the following time course; bDmin, 10s, 30s, 90s, 270s, 810s, 2430s, 7290s and Dmax. The injection and sample thawing was performed in a robotic module equipped with an universal robotics UR5 arm designed to present samples in series to a LEAP CTC-PAL system in a refrigerated cabinet. Samples were thawed and injected into the injection port following digestion. The resulting proteolytic peptide mixtures were desalted and separated on a C18 analytical column and detected by mass spectrometry. HDX-MS data treatments were performed with HDX-Workbench software. [4-5]

Results

Computational studies permitted the generation of a canonical model of the IL-23 receptor complex. To validate such model and to gain deeper insights of the interactions (dynamic view of conformational changes) between a protein and its corresponding protein partner in solution, HDX-MS experiments were performed. Thus, the proteins were exposed to D2O buffer for different periods of time to measure the deuterium incorporation in IL-23, IL-23/IL-23R, IL-23/IL-12R β 1 and IL-23/IL-23R/IL12R β 1 complexes. The HDX analysis of the Apo protein revealed that both subunits (p19 and p40) were dynamic in solution. Helices A and C and the A-B, B-C and C-D loops exhibited the highest deuterium uptake while helices B and D were identified as the less flexible regions (solvent protected) in this subunit. While some parts of domains D1-D3 in p40 exhibited faster and higher deuterium uptake, other regions in these domains showed less deuterium exchange, consequently suggesting more rigid structures. The HDX data of the IL-23/IL-23R complex indicated protection across helices A-D which allowed confirmation of the contact interactions of IL-23p19 when mixed with IL-23R. In a similar way, differential HDX analysis of IL-23 in the presence and absence of IL-12R β 1 revealed protection from exchange in mainly some parts of domains D1 and D3 of IL-23p40. The interactions of the cytokine with its partner receptors were confirmed with the HDX results in the full trimeric complex IL-23/IL-23R/IL-12R β 1.

Our results revealed that HDX-MS was also able to identify the binding epitope of phage display peptides previously identified as inhibitors of IL-23. The characterization of the interactions of

these chemotypes by HDX-MS will also be presented.

Conclusions

In conclusion, Our HDX-MS and computational analyses revealed valuable information about the interactions of IL-23 and its cognate receptors IL-23R and IL-12R β 1. These studies provide new insights into the quaternary structure of this protein-protein interaction in solution and allowed us to discover novel anti-IL-23 peptides from phage display approach. The binding activity of these peptides was confirmed by biophysical techniques meanwhile the binding epitope against IL-23 was measured by HDX-MS.

To our knowledge, this is the first structural evidence of the interactions of the full trimeric complex IL-23/IL-23R/IL-12R β 1. In addition, we confirmed by HDX-MS the characterization of the binding mode of macrocyclic small molecule to IL-23 for the first time. Indeed, this new chemotype can be used as foundation for discovery and characterization of small molecules targeting this cytokine and its receptors.

Novel Aspect

HDX-MS and computational approaches to probe the binding interactions of IL-23 with its receptors and phage display peptides.

References

1. Campa M. *et al.* *DermatolTher (Heidelb)*, 6, 1-12 (2016).
2. Espada A. *et al.* *J. Med. Chem.*, 59, 2255-2260, (2016).
3. Ting JP. *et al.* *Plos One.*, 13, 1-18, (2018).
4. Pascal B. *et al.* *J. Am. Soc. Mass Spectrom.*, 23, 1512-1521 (2012).
5. Cummins D. *et al.* *Anal. Chem.*, 88, 6607-6614, (2016).

ThOr.02-3 **MAPPING CONFORMATIONAL DYNAMICS TO BINDING MODES IN TEM-1 B-LACTAMASE USING MILLISECOND SCALE ESI-MS**

Ruth Knox - Derek Wilson

York University, Chemistry, Toronto

Keywords: catalysis, antimicrobial resistance, enzymes, HDX, beta-lactamase

Introduction

Current studies contribute to the work on antibiotic resistance by studying the enzymatic digestion of antibiotic drug families and mapping drug potency based on dynamic differences in drug binding modes. The induced fit model of enzyme function highlights the critical role of dynamics in substrate binding. However, there has been little evidence to support the notion that individual dynamic modes are required for catalytic processes after initial binding.

Methods

Tandem mass spectrometry (MS/MS) is utilized to identify peptides after proteolytic digest within a chip. Without MS/MS, column chromatography would be used to identify the eluted peptides, leaving the experiment vulnerable to back exchange and unreliable dynamic information. Experiments are performed on a Synapt G2 Time-of-Flight instrument with Ion Mobility Separation (IMS). Microfluidic devices are constructed in house using a VersaLaser and machine drill.

Results

By enabling deuterium uptake to be monitored on a millisecond timescale, we greatly enhance the manner that catalytic activity can be observed under native conditions(1). We use HDX-MS to identify dynamic modes linked to catalytic processes in the antibiotic resistant enzyme TEM-1

β -lactamase(2-4). Using a very active substrate (ampicillin), a poorly deacylated substrate (cephalexin) and an off-pathway inhibitor (clavulanate), we isolate dynamics that are linked to 3 processes: substrate binding; lactam ring hydrolysis; and, deacylation. These separations translate to leads for dynamics-targeted antibiotic development, such as at the allosteric binding site at the C-terminal helix, H11 (5,6) and the link between known mutational sites and increased conformational dynamics upon drug binding(7-9). We propose that an inner β -strand, residues 245-257, is directly responsible for increased resistance to cephalosporin drugs(10) and should be conserved through a co-inhibitor to increase the potency of current on the market drugs.

Conclusions

Our results identify dynamic modes that are uniquely linked to specific processes in the catalytic (or inhibitory) mechanism. In each case, these links could be rationalized in terms of what was known about the role of particular residues or structural elements within the identified region, supporting the reliability of the approach as a means to isolate catalytically important dynamics and to reveal new submolecular targets for drug design.

Novel Aspect

Microfluidics HDX on the millisecond timescale is unlike traditional HDX that utilizes chromatography and operates on a minute to hours timescale.

References

1. Li Z., Song F., Zhuang Z., Dunaway-Mariano D., Anderson K. S. *Anal. Biochem.* **394**, 209–216 (2009).
2. Bratulic S., Gerber F., Wagner, A. *Proc. Natl. Acad. Sci.* **112**, 201510071 (2015).
3. Ehmann D. E. et al. *Proc. Natl. Acad. Sci. U. S. Am. Early Ed.* 1–6, 6 pp. (2012).
4. De Wals P.-Y., Doucet N., Pelletier, J. N. *Protein Sci.* **18**, 147–160 (2009).
5. Horn J. R., Shoichet B. K. *J. Mol. Biol.* **336**, 1283–1291 (2004).
6. Chaibi E. B., Sirot D., Paul G., Labia, R. *J. Antimicrob. Chemother.* **43**, 447–458 (1999).
7. Henquell C., Chanal C., Sirot D., Labia R., Sirot, J. *Antimicrob. Agents Chemother.* **39**, 427–430 (1995).
8. Bershtein S., Goldin K., Tawfik, D. S. *J. Mol. Biol.* **379**, 1029–1044 (2008).
9. Caporale B., Franceschini N., Perilli M., Segatore B., Rossolini G. M. *Society* **48**, 3579–3582 (2004).
10. Salverda M. L. M., de Visser J. A. G. M., Barlow, M. *FEMS Microbiol. Rev.* **34**, 1015–1036 (2010).

ThOr.02-4 USING MASS SPECTROMETRY TO EXAMINE THE CONFORMATIONAL LANDSCAPE OF PROTEINS – THE PERFECT TOOL TO REPORT ON SELF SOLVATION AND THE EFFECT OF CHARGE LOCATION ON STRUCTURE?

Perdita Barran⁽¹⁾ - Dale Stuchfield⁽¹⁾ - Lukasz Migas⁽¹⁾ - Rebecca Beveridge⁽¹⁾ - Jakub Ujma⁽²⁾ - Bruno Bellina⁽¹⁾

⁽¹⁾ The University of Manchester, Manchester Institute of Biotechnology, Manchester - ⁽²⁾ Waters Corp., Research, Wilmslow

Keywords: Variable temperature ion mobility mass spectrometry; disordered proteins

Introduction

Mass spectrometry (MS) coupled with electrospray ionisation (ESI) has been extensively applied to identify proteins and elucidate stoichiometry of protein complexes. Since desolvated species are affected by solvent conditions such as pH, buffer strength and concentration, ESI-MS is an appropriate method by which to consider the range of conformational states that proteins may occupy including natively folded, disordered, denatured and amyloid.

Methods

Rotationally averaged collision cross sections of the ionised forms of proteins, provided by the combination of mass spectrometry and ion mobility (IM-MS), are measured and used to explore conformational landscapes in the absence of solvent. Measurements are made at a range of temperatures, on two model proteins as well as on permutants of the highly disordered protein p27.

Results

We present a novel method of monitoring the effect of charge location on structures, on unfolding transitions and the stability of transient intermediates using a variable temperature ion mobility mass spectrometer (VT IM-MS)^[1] capable of measurements from 150-500K. We show with VT IM-MS, that for both ubiquitin and fully oxidised lysozyme, transient intermediate structures can be detected upon activation followed by trapping and ion mobility separation at cryogenic temperatures. We show that for permutants of p27 the location of charged groups has a dramatic influence on the structures adopted in the gas phase, and how this compares with solution phase measurements and predictions. We have obtained data providing direct evidence that restructuring of the compact states of ubiquitin and lysozyme proceeds through extremely diverse intermediate forms. In addition, our results indicate that intermediates are initially highly elongated and subsequently converge to compacted structures.

Conclusions

The processes of elongation and subsequent compaction indicates quasi-reversible folding transitions. The influence of charge location on protein structure is strong with implications for gas phase measurements. The presented novel data obtained on well studied “model” systems highlights the potential of utilizing this methodology in studying protein folding in the gas phase.

Novel Aspect

VT-IM-MS on model proteins; consideration of the effects of sequence on stability.

References

1. J. Ujma, K. Giles, M. Morris, P. E. Barran, *Anal. Chem.* **2016**, *88*, 9469–9478.

ThOr.02-5 **A COARSE-GRAINED COMPUTATIONAL WORKFLOW FOR THE ELUCIDATION OF BIOMOLECULAR STRUCTURE AND DYNAMICS BY ION MOBILITY SPECTROMETRY-MASS SPECTROMETRY**

Rebecca D'Esposito - Daniele Fabris

University at Albany, Chemistry and Biological Sciences, Albany, NY

Keywords: coarse-graining modeling; molecular dynamics simulations; IMS-MS; collisional cross section; global dynamics.

Introduction

The interpretation of ion mobility spectrometry-mass spectrometry (IMS-MS) data from large biomolecules relies on molecular dynamics (MD) simulations to model their putative structures. With the goal of enabling millisecond simulations approaching typical IMS-MS timescales, we explored possible coarse-grained (CG) schemes to replace overly expensive all-atom frameworks, while still meeting the level of detail necessary for proper structural elucidation.

Methods

Selected nucleic acid structures were modeled by using 2-, 5-, and 7-bead representations, and then submitted to MD simulations in GROMACS by using the Martini force-field. The collisional cross section (CCS) algorithms contained in MOBCAL were reparametrized according to a previously reported procedure to accommodate the pseudoatoms included in our CG

frameworks, thus enabling the seamless calculation of theoretical CCS values directly from CG models [1].

Results

This report evaluated the merits of employing CG models to capture the global dynamics of large biomolecules, which involve the conformational rearrangement of entire domains completed in the ms range typical of IMS-MS determinations. The advantages of this approach were demonstrated by performing test MD simulations that provided a 0.009 s/ns wall time for a CG model versus 680 s/ns for the corresponding all-atom version, equivalent to 0.00025 hr versus 18 hr of computing time per 100 ns of simulation. Reparametrizing the MOBCAL algorithms provided the means for obtaining a full-fledge CG workflow to enable model building, MD simulations, and CCS calculations with no need for all-atom representations. This approach allowed us to test the possibility of employing the reparametrization process to correct for gas-phase effects on biomolecular structure, which are not adequately addressed in the leading force-fields, such as the apparent compaction of structures forced to establish self-solvating interactions in a solvent-free environment.

Conclusions

For large biomolecules, the resolution and timescale of IMS-MS analysis thwart the dissection of local dynamics involving side-chain/nucleobase displacement occurring in the ns timescale, whereas they enable the recognition of global dynamics in the ms range, which may produce distinguishable conformers. Therefore, CG approaches provide an excellent solution for completing MD simulations in the proper IMS-MS scale to enable their experimental observation.

Novel Aspect

Direct application of coarse-graining computational approaches to support the analysis of the structure/dynamics of large biomolecules by IMS-MS.

References

1. Vangaveti S. *et al. Phys. Chem. Chem. Phys.* 19, 14937-14946 (2017), DOI: 10.1039/c7cp00717e.

ThOr.03 - Natural substances and their complexes

Chairs: Pascal Gerbaux, Luisa Mattoli

ThOr.03-1 Keynote: CONTEXTUALIZING METABOLOMICS DATA TROUGH MULTI-INFORMATIVE MOLECULAR NETWORKS: METABOLITE ANNOTATION AND BIOACTIVE NATURAL PRODUCTS PRIORITIZATION STRATEGIES

Pierre-Marie Allard

School of Pharmaceutical Sciences, EPGL, University of Geneva, University of Lausanne, CMU - Rue Michel Servet, 1 1211 Geneva 4 – Switzerland

Keywords: metabolomics, pharmacognosy, contextualization approaches, molecular networks, bioinformatics

Introduction

Bioinformatic approaches are shaping multiple scientific fields and such tools have particularly impacted natural products chemistry in the latest years. The ever improving capabilities of mass spectrometry analytical platforms now allow to acquire signals of quality for the hundreds to thousands analytes constituting complex biological matrices. At the moment, such metabolomics approaches are the most appropriate to characterize natural products extracts in a

comprehensive way.

Methods

In our sense, two main challenges emerge when implementing and developing such approaches. On one side, and this defly is well identified within the community, important progresses are to be done regarding metabolite identification tools and strategies.[1] On the other hand, we propose that what we define as contextualization approaches will be crucial to fully grasp the complexity of natural products and establish the discipline of pharmacognosy in the digital era.[2]

Results, Conclusions and Novel Aspect

Here we will present some of these contextualization strategies developed in our lab. After explaining the role of central tools such as molecular networking and in-silico fragmentation approaches,[3] we will present the establishment of a bioactive natural products prioritization pipeline based on the establishment of multi-informative molecular maps.[4] We will discuss application example, evaluation methods and possible further developments.

References

1. Creek, D. J.; Dunn, W. B.; Fiehn, O.; Griffin, J. L.; Hall, R. D.; Lei, Z.; Mistrik, R.; Neumann, S.; Schymanski, E. L.; Sumner, L. W.; Trengove, R.; Wolfender, J.-L. *Metabolomics* **2014**, *8*–11.
2. Allard, P.; Bisson, J.; Azzollini, A.; Pauli, G. F.; Cordell, G. A.; Wolfender, J. *Curr. Opin. Biotechnol.* **2018**, *54*, 57–64.
3. Allard, P.-M.; Péresse, T.; Bisson, J.; Gindro, K.; Marcourt, L.; Pham, V. C.; Roussi, F.; Litaudon, M.; Wolfender, J.-L. *Anal. Chem.* **2016**, *88* (6), 3317–3323.
4. Olivon, F.; Allard, P. M.; Koval, A.; Righi, D.; Genta-Jouve, G.; Neyts, J.; Apel, C.; Pannecouque, C.; Nothias, L. F.; Cachet, X.; Marcourt, L.; Roussi, F.; Katanaev, V. L.; Touboul, D.; Wolfender, J. L.; Litaudon, M. *ACS Chem. Biol.* **2017**, *12* (10), 2644–2651.

ThOr.03-2 **NEW APPROACH FOR THE STRUCTURE ELUCIDATION OF (ACYL)POLYAMINES IN SPIDER VENOMS BY UHPLC-HR-ESI-MS/MS**

Yvonne Forster - Stefan Bienz - Laurent Bigler

University of Zurich, Department of Chemistry, Zurich

Keywords: acylpolyamine, spider venom, H/D exchange, database, in-silico fragmentation

Introduction

Acylpolyamines, often found in spider venom, are promising lead compounds to treat cancer or neurodegenerative diseases like Alzheimer's disease, epilepsy or stroke. [1,2] Polyamines follow characteristic fragmentation rules; therefore MS/MS is ideal for their structure elucidation. [3] However, low available sample amounts, co-eluting isomers and diverse nomenclature strategies in literature make the structure elucidation demanding.

Methods

The fragmentation rules for linear polyamines were used to develop a new database with in-silico fragmentation for known and unknown acylpolyamines. The spider venom was analyzed by UHPLC-HR-ESI-MS and MS/MS on a QExactive. The method was optimized for small, basic molecules. In addition, the number of acidic protons was determined in an on-line H/D exchange experiment (HDX).

Results

The database of spider venom acylpolyamines provides fast access to the related literature and

acquired MS/MS spectra. Therefore, already known acylpolyamines are identified easily and time can be used to focus on the structure elucidation of unknown lead compounds. Database entries for new acylpolyamines can easily be generated by assembling the structure from simple building blocks. The separation of co-eluting isomers could be improved. The HDX experiment clarified the presence of N-methylated amines and supports the structure elucidation in general. The optimized method allowed the identification of more than 90 polyamines from just one microgram lyophilized venom of the spider *Agelenopsis aperta*. More than half of the structures were not described before.

Conclusions

The improved analytical method in combination with the developed acylpolyamine database enables the routinely identification of linear polyamine amides even from sources of limited sample amounts like spider venoms. The in-silico fragmentation enables the identification of polyamine structures without the analysis of expensive reference material.

Novel Aspect

The database with in-silico fragmentation supports the annotation of known and the structure elucidation of new acylpolyamines based on HR-MS/MS and HDX.

References

1. Blagbrough I. S., Carrington S., Geall A. J., *J. Pharm. Sci.*, 3, 223 (1997).
2. Estrada G., Villegas E., Corzo G., *Nat. Prod. Rep.*, 24, 145 (2007).
3. Tzouros M., Chesnov S., Bigler L., Bienz S., *Eur. J. Mass Spectrom.*, 19, 57, (2013).

ThOr.03-3 GLYCOALKALOIDS IN PLANTS USING HIGH RESOLUTION MASS - NEUTRAL LOSS APPROACH

Tiziana Nardin - Roberto Larcher

Fondazione E. Mach, CTT-UCVA, San Michel a/A

Keywords: plant metabolites, Orbitrap™, liquid chromatography

Introduction

Alkaloids (alks) are widespread in nature. These nitrogen-containing organic constituents occur mainly in plants as secondary derivatives of the amino acid metabolism, or arising from amination of other substrates, which may be, for example, terpenes or steroids [1]. The biosynthesis often involves a further glycosylation of the alks (aglycones) as glycoalkaloids (glyalks) [2]. One of the best-known glyalks is the steroidal Solanine, a poisonous compound present in the Solanaceae plant family, but, to date, only a few studies on glyalks belonging to other chemical categories have been reported.

Methods

The chromatographic separation was performed using a UHPLC equipped with an SPE on-line system and a biphenyl column. The binary mobile phase was composed of 0.1% formic acid (FA) with 5 mM ammonium acetate (AAc) and MeOH/ACN 95:5 v/v with 0.1% FA and 5 mM AAc. A Full MS/AIF/NL dd-MS² experiment was performed in positive ion mode with the resolution set at 140,000 FWHM (m/z 200; 1.5 Hz) for full MS spectra, at 70,000 FWHM (3 Hz) for AIF and at 17,500 FWHM (12 Hz) for dd-MS².

Results

Neutral losses of pentose (m/z 132.0423), deoxyhexose (146.0579), and hexose (162.0528), and of all the combinations of up to four of these sugar units were considered. Full MS and AIF spectra were processed comparing the possible glyalks masses with the corresponding aglycone mass of about 100 alks already reported in literature as characteristic of the examined herbs

[3,4]. If both the glycoside (in full MS) and the aglycone (in AIF spectra) peaks were present at the same retention time, a manual interpretation of spectra was performed. The compound identification was assessed considering the matching of glycosidic experimental dd-MS² spectra with the fragmentation patterns of previously studied alks. The study allowed us to describe the profiles of glyalks belonging to the most important chemical classes (indole, piperidine, protoalkaloid, pyridine, pyrrolidine, pyrrolizidine, quinoline, steroidal, terpenoid, and tropanealks) in a selection of 21 plants.

Conclusions

This high resolution mass/neutral loss experiment permitted us to describe the presence of the glycosidic forms of a large selection of alks, permitting higher awareness of the possible risks and benefits relating to the consumption of popular plant-based medicines with new attention to both the free and bound alk forms.

Novel Aspect

To our knowledge, this is the first high-resolution mass - neutral loss approach proposed for a broad identification of plant glyalks belonging to the main chemical classes.

References

1. Aniszewski T., *Alkaloid Chemistry, Biological Significance, Application and Ecological role*, 7-10 (2007).
2. Itkin M., Rogachev I., Alkan N., Rosenberg T., Malitsky S., Masini L., Meir S., Iijima Y., Aoki K., de Vos R., Prusky D., Burdman S., Beekwilder J., Aharoni A., *The Plant Cell*, 23, 4507-4525 (2011).
3. Nardin T., Piasentier E., Barnaba C., Larcher R., *J Mass Spectrom*, 51, 729-741 (2016).
4. Nardin T., Piasentier E., Barnaba C., Larcher R., *Drug Test Anal*, 10, 423-448 (2017).

ThOr.03-4 **STRUCTURAL INVESTIGATION OF NATURAL ORGANIC MATTER USING DIRECT FRACTIONATION AND SELECTIVE ISOTOPIC LABELING COUPLED TO FTICR MS**

Alexander Zhrebker⁽¹⁾ - **Yury Kostyukevich**⁽¹⁾ - **Evgeny Shirshin**⁽²⁾ - **Oleg Kharybin**⁽¹⁾ - **Alexey Kononikhin**⁽³⁾ - **Irina V Perminova**⁽⁴⁾ - **Evgeny Nikolaev**⁽¹⁾

⁽¹⁾ Skolkovo Institute of Science and Technology, Translational medicine, Moscow - ⁽²⁾ Lomonosov Moscow State University, Physics, Moscow - ⁽³⁾ Institute for Energy Problems of Chemical Physics of RAS, Mass-spectrometry, Moscow - ⁽⁴⁾ Lomonosov Moscow State University, Chemistry, Moscow

Keywords: FT ICR, natural organic matter, isotopic exchange, fractionation, structure

Introduction

FTICR MS is an exclusive method, which enables exploration of extremely complex natural organic matter (NOM). Routine analysis resolves thousands of C_cH_hO_oN_nS_s molecular compositions in NOM sample [1]. However, FTICR MS is tolerant to structural isomers. Here we show determination of specific structural sites in individual NOM molecules and isomers in different samples using selective labeling reactions and direct fractionation in combination with FTICR MS.

Methods

NOM samples were sealed with NaOD or DCl for skeletal HDX[2]. Deuteromethylation was performed by CD₃OD/SOCl₂ reaction of NOM. All samples were desalted [3] and fractionated using SPE on ppl cartridges operated at different pH. HDX of mobile protons was performed by incubation in D₂O. All samples were analysed using 7T FT MS Bruker Apex Ultra with harmonized cell (Bruker Daltonics). The FTICR MS data were processed using the lab-made "Transhumus" software.

Results

Catalytic HDX coupled to FTICR MS was applied to determine structural fragments of individual molecules of NOM and its synthetic analog. Application of HDX showed isomeric lignin-like components, which differ by aromatic ring substitution pattern. Combination of skeletal HDX provides the reliable identification of isomers in various NOM samples isolated even from similar sources. In case of models obtained by oxidative condensation of phenols, the exact structure of components was suggested.

Deuteromethylation enabled determination of COOH-groups carried by individual molecules in NOM. This allowed for mapping compounds with close elemental compositions but different protolytic properties on conventional Van Krevelen diagram. The number of COOH-groups matched model structures suggested for major NOM components. Based on this result we performed direct fractionation of NOM on hydrophobic resin performed at different pH. We observe distinct shift of molecular composition from COOH-depleted molecules to polycarboxylic compounds.

Conclusions

We found that shared molecular compositions identified in different NOM samples correspond to structural isomers. By visualization of labeling results, we discovered relationships between compartments of NOM – oxidation, decarboxylation, condensation, etc., including transformation of CHOS compounds. Due to distinct separation of carboxylic acids according to molecular compositions, the deeper fractionation of hydrophobic components of NOM was performed.

Novel Aspect

The developed approaches enabled for the first time mapping of NOM components with respect to their structural pattern, and prediction of their chemical properties

References

1. Hertkorn, N., Ruecker, C., Meringer, M. et al. *Analyt. Bioanal. Chem.* 389, 1311-1327 (2007)
2. Zherebker A., Airapetyan D., Konstantinov A.I. et al. *Analyst.* 140, 4708-4719 (2015)
3. Dittmar T., Koch B., Hertkorn N. et al. *Limnol Ocean Methods.* 6, 230–235 (2008)

ThOr.03-5 EVALUATION OF COSMETIC INTEREST OF THE YELLOW-LEGGED HORNET VENOM (*VESPA VELUTINA NIGRITHORAX*)

Thao Nhi Le⁽¹⁾ - David Da Silva⁽¹⁾ - Cyril Colas⁽²⁾ - Eric Darrouzet⁽³⁾ - Patrick Baril⁽²⁾ - Lucie Petit Leseurre⁽⁴⁾ - Benoît Maunit⁽¹⁾

⁽¹⁾ Institut de Chimie Organique et Analytique (ICOA), UMR 7311, Université d'Orléans-CNRS, Orléans - ⁽²⁾ Centre de Biophysique Moléculaire (CBM), UPR 4301, Université d'Orléans-CNRS, Orléans - ⁽³⁾ Institut de Recherche sur la Biologie de l'Insecte (IRBI), UMR CNRS 7261, Université de Tours, Tours - ⁽⁴⁾ CHIMEX, Groupe L'Oréal, Le Thillay

Keywords: venom, antioxidant, peptidome, LC-MS/MS

Introduction

Our study aims to explore novel bioactive molecules from the Yellow-legged hornet venom for preventing or reducing skin aging in the cosmetic.

In this study two challenges are set: i) evaluate the antioxidant potency venom and then identify the bioactive structure responsible for it ii) establish and compare the peptidome expression pattern between the winter and summer venoms to assess if seasonal sampling impact qualitatively and quantitatively venom compounds.

Methods

We develop a technical investigation platform that includes DPPH (for radical scavenging capacity), FRAP (for ferric reducing antioxidant power) chemical assays as well as ROS (reactive oxygen species) production in HaCaT keratinocyte cells in order to determine the antioxidant capacity of hornet venom. To identify the active molecules responsible for this remarkable property, the crude venom was fractionated by HPLC and then identified by UHPLC-ESI-qToF MS/MS.

Results

Our results revealed a significant antioxidant activity of the venom in three activity assays. Nine fractions were obtained from HPLC and again evaluated. Only one fraction revealed an antioxidant activity. Although the fraction is rich, the identification of the compound responsible for the activity could be achieved and validated by an antioxidant test directly on Thin Layer Chromatography (TLC). From the analysis in mass spectrometry, we identified one major compound responsible for almost the whole antioxidant activity of the crude venom.

In parallel, the venom peptidomes of winter and summer hornet workers were investigated with UHPLC-qToF-HRMS or MS/MS and their differences were evaluated by statistical analysis (PCA, PLS-DA).

Conclusions

We have successfully identified the compound responsible for antioxidant activity in Yellow-legged hornet venom and compared the profile variation of venom samples issued from different periods.

Novel Aspect

The other cosmetic applications will be carried out to obtain a complete view of Yellow-legged hornet potency.

ThOr.04 - Environmental MS

Chairs: Encarnation Moyano, Cecilia Bergamini

ThOr.04-1 Keynote: EVALUATION OF COSMETIC INTEREST OF THE YELLOW-LEGGED HORNET VENOM (VESPA VELUTINA NIGRITHORAX). EMERGING CONTAMINANTS: STATE OF THE ART AND NEW DISCOVERIES

Susan D. Richardson

Department of Chemistry & Biochemistry, University of South Carolina, Columbia, SC 29205, USA

Environmental research continues to expand beyond traditional, regulated contaminants to emerging contaminants, such as sucralose and other artificial sweeteners, nanomaterials, perfluorinated compounds, pharmaceuticals, hormones, drinking water and swimming pool disinfection by-products (DBPs), 1,4-dioxane, sunscreens/UV filters, flame retardants, benzotriazoles, naphthenic acids, algal toxins, and new contaminants on the horizon: ionic liquids and microplastics. These are now frequently being found in water samples, including rivers, lakes, ground water, and drinking water.

Moreover, understanding their fate and transport in the environment and in wastewater/drinking water treatment is vitally important, and as such, one of the major trends continues to be in identifying their transformation products. Because environmental samples are inherently complex mixtures with trace-level contaminants, mass spectrometry has been key for their identification and measurement.

This presentation will provide an overview of the state of the science for emerging contaminants, their formation and transformation in the environment, and the modern tools used to measure them.

ThOr.04-2 DEVELOPMENT AND OPTIMIZATION OF A MEMBRANE-INLET-PHOTOIONIZATION MASS SPECTROMETER FOR REAL-TIME ANALYSIS OF (POLY)AROMATIC COMPOUNDS IN AQUATIC SYSTEMS

Christian Gehm⁽¹⁾ - **Thorsten Streibel**⁽¹⁾ - **Johannes Passig**⁽¹⁾ - **Detlef Schulz-Bull**⁽²⁾ - **Ralf Zimmermann**⁽³⁾

⁽¹⁾ Joint Mass Spectrometry Center, University of Rostock / Chair of Analytical Chemistry, Rostock - ⁽²⁾ Leibniz-Institute for Baltic Sea Research Warnemünde, Section of Marine Chemistry, Rostock - Warnemünde - ⁽³⁾ Helmholtz Zentrum München – German Research Center for Environmental Health, Cooperation Group Comprehensive Molecular Analytics, Neuherberg

Keywords: resonance enhanced multiphoton ionization, membrane introduction, mass spectrometry, (poly)aromatic hydrocarbons, marine environment

Introduction

Fast and sensitive analysis of (poly)aromatic hydrocarbons (PAH) is of high importance due to their impact on human health and the environment. However, very low concentrations of PAHs in marine environments complicate their detection and monitoring. A promising approach, regarding to direct determination of PAHs in marine systems, is resonance-enhanced multiphoton ionization (REMPI) coupled to membrane introduction mass spectrometry (MIMS).

Methods

In MIMS analytes are extracted from the water phase into the vacuum of the MS by using semipermeable membranes. No laborious enrichment techniques are needed. With REMPI ions are formed by the subsequent absorption of two photons via an excited molecular state. The unique ionization scheme of REMPI provides selective and sensitive detection of PAHs. For the analysis of water samples, a REMPI-MIMS system with external inlets for sheet or hollow fiber membranes was built.

Results

In this proof of principle study we present measurements of selected small (poly)aromatic species in different simulated (tap water, artificial sea water) and real world water samples. Both inlet designs show promising results with respect to the fast and sensitive determination of small PAHs. For a hollow fiber membrane setup, the trap-and-release technique can be utilized for direct determination of small (poly)aromatic species in water samples down to ppb range in less than 5 minutes. Here, the pervaporating compounds are trapped inside of the membrane tube for a certain time and are released by fast heating of the membrane as a sharp peak into the MS. By using a sheet membrane assembly fast response times and low limits of detection were obtained for toluene, *p*-xylene and naphthalene. Hereby, the membrane is heated indirectly by tempering the inflowing water. With increasing water temperature improvement of response time and signal heights can be achieved. For both inlet designs limits of detection can be improved by using high repetition lasers.

Conclusions

REMPI-MIMS is a powerful technique for the direct measurements of PAHs in aquatic environments. With this early-stage system concentrations down to the mid-ppb range for selected small (poly)aromatic compounds are easily accessible in minutes without any sample preparation.

Novel Aspect

Hyphenating external inlets with REMPI-MS is utilized for the first time. Additionally, we present the first measurements of real world water samples with REMPI-MIMS.

ThOr.04-3 SELECTIVE AND RAPID QUANTIFICATION OF ATMOSPHERIC MONOTERPENES

Kristýna Sovová - Pavel Pásztor - Violetta Shestivska - Anatolii Spesyvyi - Michal Lacko - Patrik Španěl

J. Heyrovský Institute of Physical Chemistry of the CAS, v. v. i., Department of Chemistry of Ions in Gaseous Phase, Prague

Keywords: SIFT-MS, SIFDT-MS, VOCs, TD, FastGC, quantification

Introduction

The main objective is to develop a new analytical method for selective and rapid quantification of individual monoterpene isomers. This is especially relevant for atmospheric chemistry because different isomers exhibit highly variable reactivity with the OH[•] radical [1] and their role in formation of pollutants may be different for individual isomers. Thus, the precise measurements of their fluxes and concentrations will refine input data for numerical models.

Methods

Soft chemical ionization technique of SIFT-MS [2] and its drift tube variant [3] have been used for absolute quantification of monoterpene isomers. Because SIFT-MS can quantify only the total concentration of monoterpene isomers, we employed more traditional separation method of fast gas chromatography, GC, [4] and thermal desorption, TD, and coupled them with *Profile 3* instrument.

Results

Separation by the fast GC was carried out on an electrically heated 5 m long non-polar MXT capillary column. Elution times of the isomeric compounds through short, rapidly heated column were characterised and optimised. Using the fast GC coupled to a *Profile 3* SIFT-MS instrument we successfully separated eight of the most abundant plant monoterpene isomers (α -pinene, β -pinene, Camphene, Myrcene, 3-carene, R-limonene, α -terpinene and γ -terpinene). A combination of TD with SIFT-MS showed the ability for partial separation of isomers achieved by a slow temperature desorption profiles in contrast to the common approach of rapid desorption from TD tubes. The characterization of several sorbents was experimentally researched and tested. The optimal conditions for the TD approach such as the amount of sorbent, sampling flow rate, or sampling volume were defined.

Conclusions

The combination of fast GC or TD with SIFT-MS represents a promising and robust analytical method for selective and rapid quantification of monoterpene isomers. The advantage of fast GC is better separation in time of a minute, can be used in-situ, but does not allow the use of eddy covariance method. TD based method is very useful for long term monitoring of atmosphere and sample storage on TD tubes.

Novel Aspect

Separation by TD profiles of the isomeric compounds trapped on sorbent was researched. Separation by fast GC using short capillary column was characterised and optimised.

References

1. K. Hens, A. Novelli, M. Martinez, J. Auld, R. Axinte, B. Bohn, *Atmospheric Chemistry and*

- Physics*, 14, 8723-8747 (2014).
2. D. Smith, P. Španěl, *Mass Spectrometry Reviews*, 24661-700 (2005).
 3. A. Spesyvyi, D. Smith, P. Spanel, *Analytical Chemistry*, 87 12151-12160 (2015).
 4. D. Materić, M. Lanza, P. Sulzer, J. Herbig, D. Bruhn, C. Turner, N. Mason, V. Gauci, *Analytical and Bioanalytical Chemistry*, 407, 7757-7763 (2015).

ThOr.04-4 **COMPREHENSIVE DEPTH PROFILING OF STATEN ISLAND SOIL BY GC AND PETROLEOMICS**

Mary Thomas⁽¹⁾ - **Emma Collinge**⁽²⁾ - **Matthias Witt**⁽³⁾ - **Christopher Vane**⁽⁴⁾ - **Mark P. Barrow**⁽²⁾

⁽¹⁾ MAS CDT, University of Warwick, Coventry - ⁽²⁾ Department of Chemistry, University of Warwick, Coventry - ⁽³⁾ Bruker Daltonik GmbH, Bruker Corporation, Bremen - ⁽⁴⁾ Organic Geochemistry, British Geological Survey, Keyworth

Keywords: Contaminated soil, petroleomics, GC, FTICR MS

Introduction

Staten Island is located in the New York/New Jersey (NY/NJ) Estuary, a region with one of the highest population densities in the United States. Surface sediment concentrations showed the NY/NJ Estuary to be among the most chemically contaminated waterways in the United States. Analysis of the extent of environmental contamination by complex petroleum-related pollutants is of rising importance.

Methods

Sediment cores were collected in a marsh region of Staten Island, NY and freeze-dried. Samples taken at intervals to 1 m depth were spiked with PAH standard mixture and subjected to solid phase extraction (SPE) before bulk analysis by GC MS. Sediment samples from 5 selected depths underwent Soxhlet extraction prior to compositional analysis by direct infusion APPI-FTICR MS and GC-APCI-FTICR MS.

Results

Saturate, aromatic and resins (SAR) GC MS analysis at each sampling depth showed a sharp increase in concentration at two depths that may correlate with the timing of major spills in the area, including the 1990 Exxon pipeline spill. Direct infusion APPI-FTICR MS measurements showed a corresponding spike in the relative contributions to spectral intensity from radical hydrocarbon as well as oxygen and sulfur-containing compound classes at these depths. However, background sampling depths were also found to have strong contributions from petroleum-related contaminants, indicating continuously poor water and soil quality in the NY/NJ Estuary resulting from regular oil spills in the region. Using GC-APCI-FTICR MS additional compound classes were observed, revealing the impact of leaks of industrial chemicals and agricultural effluent on the estuary system.

Conclusions

Several trends in the compositional profiles obtained by direct infusion APPI-FTICR MS and GC-APCI-FTICR MS were found to correlate with the bulk analysis obtained by GC-MS, including relatively strong contributions from O_xS_y[H] compound classes. Samples originally obtained for background study were found to contain significant contamination from anthropogenic sources.

Novel Aspect

Comprehensive study of petroleum-contaminated soil by GC-MS, direct infusion APPI-FTICR MS and GC-APCI-FTICR MS.

ThOr.04-5 **NEW PESTICIDE METABOLITES IN GROUNDWATER DETECTED USING HIGH RESOLUTION MASS SPECTROMETRY**

Karin Kiefer - Adrian Müller - Heinz Singer - Juliane Hollender

Eawag, Environmental Chemistry, Dübendorf

Keywords: suspect screening, high resolution mass spectrometry, pesticide metabolites, groundwater

Introduction

Groundwater contributes to 80% of drinking water supply in Switzerland [1]. However, various contaminant sources such as pesticide application in agriculture pose risks to groundwater quality. Only for few pesticide metabolites reference standards are commercially available so that groundwater quality is partly unknown. Thus, high resolution mass spectrometry combined with a screening of exact masses of potential metabolites is a possibility to detect new contaminants.

Methods

31 groundwater samples were enriched via water evaporation, and analyzed with liquid chromatography coupled with electrospray ionization to high resolution tandem mass spectrometry (LC-ESI-HRMS/MS). The HRMS data were searched for the exact masses of more than 1000 pesticide metabolites, most of them observed within the European pesticide registration process [2, 3], using Compound Discoverer (Thermo Scientific).

Results

After background and target removal 9300 compounds were annotated as suspects in positive and negative ionization mode. Despite optimization of the parameter settings in Compound Discoverer using target compounds, a high number of compounds were regrouped insufficiently across samples and noise peaks were not fully removed. Therefore, visual inspection was required. The remaining 680 compounds, which showed a signal/noise >10, were checked for plausibility based on isotope pattern, retention time and in silico fragmentation with MetFrag [4]. The pesticide producers provided reference standards for 28 metabolites originating from 16 pesticides. So far, 11 compounds have been confirmed; 4 out of the 11 compounds were already detected in a target screening study by Reemtsma et al. [5].

Conclusions

Due to the high-quality suspect list, the optimized workflow and the HRMS data the suspect screening successfully identified several new pesticide metabolites which have not been measured before in groundwater. However, despite advances in the evaluation softwares during recent years, suspect screening is still related to a large manual effort, a combination of computational tools and expert knowledge.

Novel Aspect

Groundwater samples were screened for >1000 pesticide metabolites. Water evaporation allowed to minimize losses of polar compounds during enrichment.

References

1. Tripet J., *The hydrology of Switzerland: selected aspects and results*, 79-100 (2005).
2. Latino D.A., Wicker J., Gutlein M., Schmid E., Kramer S., Fenner K., *Environmental Science: Processes & Impacts*, 19, 449-464 (2017).
3. Lewis K.A., Tzivilakis J., Warner D.J., Green A., *Human and Ecological Risk Assessment: An International Journal*, 22, 1050-1064 (2016).
4. Ruttkies C., Schymanski E.L., Wolf S., Hollender J., Neumann S., *Journal of Cheminformatics*, 8, 3 (2016).

5. Reemtsma T., Alder L., Banasiak U., *Water Research*, 47, 5535-5545 (2013).

ThOr.05 - Safe & Valuable Food

Chairs: Michel W.F. Nielen, Emanuela Gregori

ThOr.05-1 Keynote: HOW FAR CAN WE TRUST UNTARGETED HRMS-BASED EXPOSOMICS TO ASSESS FOOD VALUE AND SAFETY FOR THE CONSUMER: THE CASE OF EXPOSURE TO PESTICIDES

Laurent Debrauwer

INRA, UMR 1331 TOXALIM, TOULOUSE, France

Keywords: food safety, metabolomics, pesticides, high resolution

Introduction

The presence of chemical contaminations in the environment and food can lead to human exposure to various molecules (including pesticides), which could be responsible for adverse effects for the Human consumer. Biological monitoring allows improving the exposure assessment, but this usually concerns a limited number of compounds analyzed in a targeted way, while untargeted approaches should provide a wider and more relevant assessment of this exposure.

Methods

HRMS data were acquired on LTQ-Orbitrap XL (Thermo) and Synapt G2-SI(Waters) mass spectrometers running under both positive and negative electrospray ionization modes. Both spectrometers were fitted with UHPLC systems using generic H₂O - MeOH gradient elution on classical C18 columns. Data were treated using various software solutions including MassLynx and Progenesis from Waters, Xcalibur, MetWorks and TraceFinder from Thermo Scientific, XCMS, SIMCA-P.

Results

A non-targeted UHPLC-ESI-HRMS exposomics workflow has been developed [1] and applied to the analysis of Human cohort urine samples. After HRMS detection, data can be processed in two ways: (i) focusing on a list of suspected metabolites generated from selected contaminants (e.g. pesticides), or (ii) starting from all features for identification of discriminant metabolites (including both endogenous and xenobiotics metabolites). From the same dataset, the first approach allows untargeted exposure assessment whereas the second can highlight possible impact of these exposures on the health status.

Through few examples, this lecture will give illustrations of how untargeted HRMS allows characterisation of pesticide metabolites as well as other chemicals, representing markers which enable classifying various population groups [2]. Other complementary strategies can also be implemented in order to get more toxicological relevance, by focusing e.g. on toxic reactive metabolites formed in-vivo.

Conclusions

The combination of untargeted and suspect screening approaches, constitutes an efficient tool for exposomics studies, and could be used as an upstream multi-exposure assessment step for food and environmental exposure evaluation. The relevance and limitations of those untargeted or semi-targeted approaches for better assessment of food safety in the context of multi-exposure to contaminants will be discussed.

Novel Aspect

Untargeted and suspect screening data acquisition strategies for evidencing potentially new biomarkers of exposure to chemical contaminants

References

1. Jamin E.L., Bonvallet N., Tremblay-Franco M., Cravedi J.P., Chevrier C., Cordier S., Debrauwer L., *Analytical and Bioanalytical Chemistry*, 406, 1149 (2014).
2. Jamin E.L., Delcambre A., Martin J.F., Bonvallet N., Kesse-Guyot E., Cravedi J.P., Debrauwer L., *Proc. 65th ASMS Conference, Indianapolis, IN* (2017).

ThOr.05-2 NOVEL IDENTIFICATION OF 3-EPI-25-HYDROXYVITAMIN D₃ IN PORK AND BEEF LONGISSIMUS THORACIS STEAKS CONFIRMS THE PRESENCE OF THIS METABOLITE IN THE HUMAN FOOD CHAIN

Kirsten G. Dowling⁽¹⁾ - Aoife Hayes⁽¹⁾ - Karen Galvin⁽¹⁾ - Sarah K. Duffy⁽²⁾ - Alan K. Kelly⁽²⁾ - John V. O'Doherty⁽²⁾ - Kevin D. Cashman⁽¹⁾

⁽¹⁾ Cork Centre for Vitamin D and Nutrition Research, School of Food and Nutritional Sciences, University College Cork, Cork - ⁽²⁾ School of Agriculture and Food Science, University College Dublin, Dublin

Keywords: 3-epi-25(OH)D₃, 25-hydroxyvitamin D₃, beef, pork, LC-MS/MS

Introduction

The source of the C3 epimer (3-epi-25(OH)D₃) of 25-hydroxyvitamin D₃ (25(OH)D₃) in human serum may arise, in part, from exogenous intake [1]. While pork and beef contain 25(OH)D₃ [2,3], it is not known if pork and beef supply 3-epi-25(OH)D₃ to the human diet. Measurement of 3-epi-25(OH)D₃ in pork and beef and in the corresponding serum of domestic animals may provide clues to the epimer's origin.

Methods

A UPLC-MS/MS method in serum [4] was adapted for the analysis of 25(OH)D₃ in *Longissimus thoracis* steaks (LT). The pentafluorophenyl (PFP) column was able to separate isobaric 3-epi-25(OH)D₃ from 25(OH)D₃ in both matrices. Quantification of 3-epi-25(OH)D₃ in LT was achieved using 25(OH)D₃-d₆ as the IS and applying a correction factor to adjust for different ionization efficiencies of the epimers.

Results

The LOD for 3-epi-25(OH)D₃ in LT was 0.008 ng/g, and the RSD% of minced beef control (n=13 over 7 days) was 30%. The concentrations of 3-epi-25(OH)D₃ in LT (n=18 pork, n=19 beef) ranged from 0.016-0.670 ng/g. There was no significant difference between mean 3-epi-25(OH)D₃ concentrations in pork and beef (P>0.05), and the overall mean±SD was 0.102±0.070 ng/g. Thus, a typical 150g pork or beef steak would provide 15.0 ng 3-epi-25(OH)D₃. In serum collected from the same animals, the mean±SD 3-epi-25(OH)D₃ was higher in porcine (2.9±0.9 ng/mL) compared to bovine serum (1.9±0.7 ng/mL). 3-epi-25(OH)D₃ was positively correlated to 25(OH)D₃ in both LT (R=0.570) and serum (R=0.444). The %epimer (ie. 3-epi-25(OH)D₃/(3-epi-25(OH)D₃+25(OH)D₃)*100) was equivalent in matched LT-serum pairs (P>0.05), and the mean±SD% epimer in beef LT was 6.2±3.5%, compared to 4.0±1.8% in pork LT.

Conclusions

This is the first known evidence of 3-epi-25(OH)D₃ in the human food chain originating from pork and beef LT steaks. The % epimer in LT was the same as that in the circulation of both species, which suggests epimerisation occurs endogenously within the animals. The difference in %

epimer content in pork compared to beef highlights species-related differences in the rate of epimerisation.

Novel Aspect

The PFP column commonly used to separate 25(OH)D₃ epimers in serum was used to identify 3-epi-25(OH)D₃ in pork and beef LT steaks.

References

1. Engelman C.D., Bo R., Zuelsdorff M., Steltenpohl T., *et al.*, *Clin Nutr*, 33, 412-425 (2014).
2. Duffy S. K., O'Doherty J. V., Rajauria G., Clarke L. C., *et al.*, *Meat Sci*, 134, 103-110 (2017).
3. Duffy S., Kelly A., Rajauria G., Jakobsen J., *et al.*, *Meat Sci* (1st revision under review).
4. Cashman K. D., Kinsella M., Walton J., Flynn A., *et al.*, *J Nutr* 2014, 144, 1050-1057 (2014).

ThOr.05-3 ARTIFICIAL INTELLIGENCE ENABLES THE DETECTION OF PATHOGENS IN FOOD AND OTHER COMPLEX BIOLOGICAL MATRICES BY MASS SPECTROMETRIC ANALYSIS OF RNA MODIFICATIONS

Botros Toro⁽¹⁾ - **Reza Nemati**⁽²⁾ - **Mehraveh Salehi**⁽³⁾ - **Will McIntyre**⁽²⁾ - **Colin Aldrich**⁽²⁾ - **Daniele Fabris**⁽⁴⁾

⁽¹⁾ University at Albany, Biological Sciences, Albany, NY - ⁽²⁾ University at Albany, Chemistry, Albany, NY - ⁽³⁾ Yale University, Electrical Engineering, New Haven, CT - ⁽⁴⁾ University at Albany, Chemistry and Biological Sciences, Albany, NY

Keywords: machine learning; pathogen detection; biomarker panels; food safety; medical diagnostics.

Introduction

We have recently demonstrated that a variety of microorganisms and viruses can be accurately differentiated according to their specific RNA post-transcriptional modifications (PTMs) observed by MS analysis [1]. We have now employed machine learning to establish a classification and identification platform that capitalizes on comprehensive PTM profiles to enable pathogen detection not only in cell cultures, but also in actual biological samples and materials.

Methods

Comprehensive PTM profiles were obtained from total RNA extracts treated with exonucleases to generate mononucleotide mixtures, which were analyzed by direct infusion nanospray MS on either an Orbitrap Velos (Thermo) or a Synapt G2 (Waters) instrument [2]. A non-redundant database of RNA modifications created in house was employed to perform data interpretation and generate final datasets used for training machine-learning algorithms written in Python.

Results

Machine learning constitutes an ideal approach for establishing a classification/identification platform in which the unique differentiating feature consists of an entire panel of biomarkers, including identities and corresponding abundances, rather than an individual biomarker. In our case, each unique panel consists of covalent ribonucleotide variants (i.e., PTMs) identified according to their characteristic masses, fragmentation patterns, and arrival time distributions. In this report, we show that the PTM profile of a typical microorganism –comprising on average ~30 PTMs– can be unambiguously retraced to its source. We demonstrate the ability of either supervised or unsupervised algorithms to properly classify labeled samples and recognize unlabeled (“unknown”) ones. Among the different systems included in the study, our classification/identification platform correctly differentiated milk contaminated by *S. typhimurium* from uncontaminated samples, as well as mosquitoes infected by Zika virus from healthy controls.

Conclusions

We demonstrated that different machine learning algorithms enable the correct identification of bacterial/viral pathogens in samples containing high backgrounds of foreign RNA, which were afforded, for example, by milk somatic cells or mosquito tissues. The ability to identify pathogens in their typical biological matrices and without analytical separations has an enormous potential for a broad range of food safety and medical diagnostic applications.

Novel Aspect

Demonstration of the ability of artificial intelligence to support pathogen identification through the analysis of large biomarker panels.

References

1. McIntyre, W., et al. *Nucl. Acids Res.* (2018) in press DOI: 10.1093/nar/gky029.
2. Rose, R.E., et al. *Mol. Cell. Proteomics* 15, 932-44 (2016) DOI: 10.1074/mcp.M115.054718.

ThOr.05-4 **HYDROGEN/DEUTERIUM EXCHANGE IS AN EFFECTIVE METHOD FOR THE STRUCTURAL RESOLUTION OF ISOMERIC PROANTHOCYANIDINS IN WINE AND OTHER DIETARY SOURCES**

Vakare Merkyte⁽¹⁾ - **Edoardo Longo**⁽¹⁾ - **Fabrizio Rossetti**⁽²⁾ - **Emanuele Boselli**⁽¹⁾

⁽¹⁾ Free University of Bozen – Bolzano, Faculty of Science and Technology, Bolzano - ⁽²⁾ Università Politecnica delle Marche, Department of Agricultural, Food and Environmental Sciences, Ancona

Keywords: isotopic exchange; high resolution mass spectrometry; proanthocyanidins, wine

Introduction

Proanthocyanidins (PAC) show high degree of isomerism, which makes it difficult to confirm the identification of new congeners. A novel class of PAC with a macro-cyclic structure has been recently found in wine [1, 2]. Hydrogen/deuterium isotopic exchange (HDX) combined to high resolution MS is discussed for the disambiguation of such cyclic B-type class from other isomeric PAC present in wine and other significant dietary sources.

Methods

HPLC-HDX-HRMS and HPLC-HRMS/MS were applied to an SPE purified and concentrated Lagrein (a red wine from South Tyrol, Italy), a peanut skin extract and a cranberry extract. For SPE purification a published methodology was followed [3]. The instrumental method was adapted from a published report [1].

Results

A series of 4-, 5- and 6-mers PAC were tentatively identified by HPLC-HRMS/MS in wine and (in lower extent) in cranberries. Their MS/MS fragmentations and retention time profiles showed remarkable differences with non-cyclic B-type analogues (already known and studied): a) beside the expected -2H difference, much lower degree of MS/MS fragmentations and very different retention time distributions distinguished the cyclic B-type PAC from the more abundant non-cyclic B-type PAC; b) tentative cyclic analogues eluted faster and showed fewer (often only one) isomeric peaks. However, cyclic PAC could not theoretically be distinguished from the isomeric non-cyclic A-type PAC. HPLC-HDX-HRMS was applied with this aim: the analysis of a peanut skin extract provided evidences that upon HDX the isomeric A-type PAC exchange one proton less than the cyclic B-types. By means of this approach, several cyclic B-type 4-, 5- and 6-mer PAC were confirmed in wine, whereas mostly non-cyclic A type were confirmed in peanut skin.

Conclusions

Isotopic exchange HPLC-HRMS proved to be a powerful tool, complementary to NMR and other

techniques, to investigate the macromolecular structure of proanthocyanidins from several plant-derived foods. In fact, a whole series of novel cyclic B-type 4-, 5- and 6-mer proanthocyanidins could be distinguished from their isomeric non-cyclic A-type analogues by isotopic exchange mass spectrometry only.

Novel Aspect

A mass spectrometric method based on hydrogen/deuterium exchange for the resolution of novel cyclic isomeric proanthocyanidins in wine and plant-derived foods is presented.

References

1. Longo E., Rossetti F., *et al.*, *Journal of the American Society for Mass Spectrometry*, 29(4), 663 (2018).
2. Jouin A., Rossetti F., Teissède P. *et al.*, *10th IVAS Symposium, Salamanca (Spain)*, 160 (2017).
3. Jeffery D., Mercurio M.D., Herderich J. *et al.*, *Journal of Agricultural and Food Chemistry*, 56(8), 2571 (2008).

ThOr.05-5 MALDI-TOF MASS SPECTROMETRY IDENTIFICATION AND DETECTION OF RELEVANT PATHOGENS IN BEEF CATTLE

Gisele Bacanelli⁽¹⁾ - **Cynthia Mantovani**⁽¹⁾ - **Anna Letícia Louzan**⁽¹⁾ - **Taynara Pasquatti**⁽²⁾ - **Daniele Bier**⁽²⁾ - **Gracia Rosinha**⁽³⁾ - **Flabio Araujo**⁽³⁾ - **Newton Verbisck**⁽³⁾

⁽¹⁾ Federal University of Mato Grosso do Sul, UFMS, Campo Grande - ⁽²⁾ Dom Bosco Catholic University, UCDB, Campo Grande -

⁽³⁾ Embrapa Beef Cattle, Animal Health, Campo Grande

Keywords: MALDI-TOF, Mycobacterium, Brucella, Salmonella, Beef Cattle.

Introduction

Mycobacterium bovis, *Brucella abortus* and *Salmonella* species account for significant economic losses in cattle production worldwide, besides being public health threats. We have established MALDI-TOF (Matrix Assisted Laser Desorption Ionization – Time-of-Flight) mass spectrometry with modified protein extraction methods and improved the presently available reference spectra libraries, in order to effectively detect those pathogens from bovine samples.

Methods

Bacterial isolation from bovine tissue and lesions followed microbiological standard methods. Heat inactivated mycobacteria from solid media were mechanically disrupted as designed herein. *Brucella* and *Salmonella* were grown in liquid media and processed as described [1]. Mass spectra acquired with alpha-cyano-4-hydroxycinnamic acid on Autoflex III Smartbeam underwent microorganism identification carried out on MALDI Biotyper software (Bruker Daltonik).

Results

Analysis of mycobacterial isolates revealed that sixty-three classified as *Mycobacterium bovis*, two as *Gordonia sputi* and one as *M. nonchromogenicum*. These results were confirmed by Polymerase Chain Reaction (PCR) specific for *M. bovis* and 16S rDNA sequencing (100% concordance).

Brucella abortus, *B. suis*, *B. ovis* and *B. canis* mass spectra profiles, generated as reference strains, were clearly different and allowed field samples identification. Amongst twelve isolates with gender specific PCR positive, only two were confirmed as *Brucella* by MALDI-TOF, indicating possible contaminations during cultivation. We also tested if *Brucella* fingerprint could be detected directly from tissue lesion homogenates, lacking any microbiological culture, which was successfully demonstrated.

Ninety-six isolates from bovine carcasses were biochemically tested for *Salmonella* and species were identified by MALDI Biotyper analysis. Four isolates showing inconclusive results in biochemistry actually belonged to *Citrobacter* and *Proteus* genera [2].

Conclusions

The cell processing method established here provides reliable *Mycobacterium* identification at species level, with scores that safely allow *M. bovis* and *M. tuberculosis* distinction of isolates from solid media culture. We have discriminated *Brucella* species and could detect this pathogen directly from bovine tissue lesions, skipping bacterial cultivation step. In addition, MALDI-TOF enabled *Salmonella* identification at species level with serovar indicative.

Novel Aspect

Distinction of Mycobacterium Tuberculosis Complex members at species level and detection of *Brucella* directly from infected tissue by MALDI-TOF mass spectrometry.

References

1. Freiwald A. & Sauer S., *Nature Protocols*, 4, 732 (2009).
2. Bier D., Tutija J., Pasquatti T., Oliveira T., Araujo F., Verbisck N., *Pesquisa Veterinaria Brasileira*, 37, 1373 (2017).

ThOr.06 – Lipidomics

Chairs: Gavin Reid, Donatella Caruso

ThOr.06-1 **Keynote:** LIPIDOMICS ANALYSIS OF CELLS AND TISSUES IDENTIFIES THERAPEUTIC TARGETS

Michael Wakelam⁽¹⁾ - **An Nguyen**⁽¹⁾ - **Qifeng Zhang**⁽¹⁾ - **Roberto Solari**⁽²⁾

⁽¹⁾ Babraham Institute, Signalling, Cambridge, United Kingdom - ⁽²⁾ Imperial College, Imperial College, London, United Kingdom

Keywords: Mass spectrometry, lipidomics, pathway analysis, colorectal cancer, rhinovirus

Introduction

Theoretically mammalian cells can contain many thousands of individual lipid molecular species, whilst it is unlikely that all are present in a single cell, lipidomic experiments have demonstrated the presence of more than a thousand species. The integrated regulation of changes in lipid species that occur that can regulate cell functions, including signaling and metabolism, highlights the need for bioinformatics analysis to fully interpret lipidomics data.

Methods

We have adopted two experimental systems: lipids extracted from colorectal tumour tissue and rhinovirus-infected human bronchial epithelial cells, were identified by LC-MS/MS. Pathway analysis of lipid metabolising enzymes coupled to network optimising Prize-collecting Steiner tree problem methodology was adopted to identify key enzymatic changes which could be considered as therapeutic targets.

Results

Lipidomic analysis that compared colorectal tumour tissue of defined stage with matched normal tissue showed changes in more than 700 lipid species, this also demonstrated changes in acyl chain saturation and chain length accompanying tumour progression. In the rhinovirus-

infected cells there was temporal-associated changes in a number of lipids including PIP, cardiolipin, PC and DG, there were also clear changes in the acylchains of ofneutral and phospholipids.

Pathway analysis highlighted key changes in sphingolipid metabolising enzymes in both systems including sphingosine kinase and ceramide synthase activities. Further, adopting the network optimising Prize-collecting Steiner tree problem methodology, we were able to identify fatty acid synthases, acyltransferases and ceramidases as potential therapeutic targets. The importance of a number of these activities in viral infection was confirmed by observing inhibitory effects of selective inhibitors upon infection of human bronchial epithelial cells by rhinoviruses.

Conclusions

The lipidomics-determined changes identify potential novel therapeutic targets, they additionally suggest roles for remodelling of membrane lipid acyl chains thus changes in membrane fluidity, protein binding or curvature may have important roles in infection and disease progression.

Novel Aspect

Lipidomics with appropriate bioinformatics analysis can point to key enzyme activities that are responsible for disease-associated change. This can therefore point to novel therapeutic targets that can be subsequently explored.

ThOr.06-2 **NON-SUBJECTIVE CHARACTERIZATION OF MEDULLOBLASTOMA SUBGROUP AFFILIATION WITH DESORPTION ELECTROSPRAY IONIZATION MASS SPECTROMETRY**

Arash Zarrine-Afsar⁽¹⁾ - ***Michael Woolman***⁽²⁾ - ***Alessandra Tata***⁽¹⁾ - ***Isabelle Ferry***⁽³⁾ - ***Claudia Kuzan-Fischer***⁽³⁾ - ***Megan Wu***⁽³⁾ - ***Sunit Das***⁽³⁾ - ***Michael Taylor***⁽³⁾ - ***James Rutka***⁽³⁾ - ***Howard Ginsberg***⁽⁴⁾

⁽¹⁾ University of Toronto, University Health Network, TECHNA Institute, Toronto - ⁽²⁾ University of Toronto, Medical Biophysics, Toronto - ⁽³⁾ Hospital for Sick Children, University of Toronto, Toronto - ⁽⁴⁾ St Michael's Hospital, University of Toronto, Toronto

Keywords: Desorption electrospray ionization mass spectrometry, medulloblastoma, lipid profiling, rapid pathology, tissue smears

Introduction

Medulloblastoma (MB) is a heterogeneous pediatric brain tumour comprised of 4 molecular subgroups that appear morphologically identical, under the microscope for conventional intraoperative, rapid pathologic assessment [1-3]. While each MB subgroup presents its own survival rate, unique response to chemotherapy and the necessary aggressiveness of resection, they remain indistinguishable intraoperatively.

Methods

Desorption electrospray ionization mass spectrometry is an effective technique for tumour type identification based on profiling of tissue lipid and small metabolites [4,5]. Here we subjected the rapidly prepared smears of murine models of medulloblastoma tumours to DESI-MS analysis to harness both the specificity and the sensitivity of this technique in determining MB subgroup affiliation in less than 1 min of sample preparation and analysis time.

Results

Mouse cell-line derived xenograft tumours of MB Group 3 and SHH subgroups were subjected to preparation of smear samples (30 seconds) and 100 seconds of DESI-MS data collection for the preparation of multivariate statistical model that can distinguish between the two subgroups using. The model was comprised of 400 data points (200 per subgroup) comprised of 5 independent cell lines and 10 independent tumours, and showed clear separation of data points

belonging to each subgroup. The model was validated using 20% leave out test. The predictive ability of the model was then assessed using blind tests. Smear spreads from two independent tumours belonging to each MB subgroup were prepared in 30 seconds and subjected to additional 10 seconds of DESI-MS data collection and real time analysis against the multivariate model using binned spectral data. A success rate of greater than 95% was achieved using the combined platform of rapid DESI-MS analysis of smears and online multivariate recognition and model comparison platforms.

Conclusions

Rapid determination of medulloblastoma subgroups based on DESI-MS lipid profiles in less than 1 min of smear sample preparation and analysis time is shown to be possible. This observation allows prognostic information regarding surgical margin extent tailored to the specific requirement of each MB subgroup to be fully utilized within the time constraints of surgery, further highlights specificity of lipid profiling with DESI-MS.

Novel Aspect

Identification, within 1 minute, of morphologically identical subgroups of same tumour type with DESI-MS analysis of tissue smears based on lipid profiling is possible.

References

1. Northcott, P. A. *et al.* Medulloblastoma comprises four distinct molecular variants. *J Clin Oncol* 29, 1408-1414 (2011).
2. Ostrom, Q. T. *et al.* CBTRUS Statistical Report: Primary Brain and Central Nervous System Tumors Diagnosed in the United States in 2008-2012. *Neuro Oncol* 17 Suppl 4, iv1-iv62, (2015).
3. Ostrom, Q. T. *et al.* Alex's Lemonade Stand Foundation Infant and Childhood Primary Brain and Central Nervous System Tumors Diagnosed in the United States in 2007-2011. *Neuro Oncol* 16 Suppl 10, x1-x36 (2015).
4. Zhang, J., Yu, W., Suliburk, J. & Eberlin, L. S. Will Ambient Ionization Mass Spectrometry Become an Integral Technology in the Operating Room of the Future? *Clin Chem* (2016).
5. Ifa, D. R. & Eberlin, L. S. Ambient Ionization Mass Spectrometry for Cancer Diagnosis and Surgical Margin Evaluation. *Clin Chem* 62, 111-123 (2016).

ThOr.06-3 EXPANDING THE LIPIDOME BY DATA-INDEPENDENT OZONE-INDUCED DISSOCIATION TO PROFILE COMPLEX LIPID ISOMERS

David L. Marshall - Berwyck Poad - Reuben Young - Stephen Blanksby

Queensland University of Technology, Institute for Future Environments, Brisbane

Keywords: lipidomics, ozone-induced dissociation, data-independent analysis, plasma

Introduction

Human plasma contains over 1500 different lipids at the sum composition level.^[1] Yet the number of unique lipids present is surely higher, once isomeric contributions from double bond location and fatty acyl regiochemistry are considered. To close this gap, ozone-induced dissociation (OzID) is incorporated into data-independent lipidomics workflows, thus enabling complex extracts to be profiled for unsaturation and regioisomerism by construction of 2D maps.

Methods

Ozone-induced dissociation was implemented on a Thermo Fisher Scientific LTQ Orbitrap Elite high-resolution mass spectrometer, equipped with an ESI source operating in positive ion mode. The instrument is modified such that nitrogen is replaced by ozone (~15% in oxygen) as the HCD

collision gas. Lipid standards and complex extracts were introduced into the ESI source by direct infusion (Advion TriVersa Nanomate).

Results

OzID performance was benchmarked using glycerolipid standards containing unsaturated fatty acyl chains. Isolation of lipid ions in the linear ion-trap, followed by interaction with ozone in the HCD cell yielded diagnostic ions for each double bond position and *sn* position in the resulting MS/MS spectrum. The increased ozone concentration in the HCD cell (compared to a single stage ion-trap instrument^[2]) significantly improved ozonolysis efficiency. This advance translates into increased lipidome coverage, and improvements in duty cycle for data-independent MS/MS analysis using shotgun workflows. Grouping all precursor ions with a common OzID neutral loss enables straightforward classification of the lipidome by unsaturation position (with respect to the methyl terminus). Glycerolipids exhibiting common HCD/OzID product ions were classified according to the fatty acyl carried at the *sn*-1 position. 2-dimensional maps provide a straightforward visualization of the data.

Conclusions

Resolving ambiguity in the molecular structure of lipids is the next frontier in contemporary lipidomics.^[3] OzID provides access to this detailed level of molecular characterization. Global profiling of fatty acyl unsaturation within a biological extract provides a snap-shot of enzymatic biosynthesis and could potentially provide indicators of health and disease.

Novel Aspect

First demonstration of a data-independent strategy combined with ozone-induced dissociation for comprehensive profiling of unsaturation in complex lipid samples.

References

1. Bowden J., Heckert A., Ulmer C., Jones C., Koelmel J *et al.*, *J. Lipid Res.*, DOI: 10.1194/jlr.M079012 (2017).
2. Thomas M., Mitchell T., Harman D., Deeley J., Nealon J., Blanksby S., *Anal. Chem.*, 80, 303 (2008).
3. Zacek P., Bukowski M., Rosenberger T., Picklo M., *J. Lipid Res.*, 57, 2225 (2016).

ThOr.06-4 **MODELLING CANCER LIPOGENESIS USING LA-REIMS METABOLIC FLUX ANALYSIS IN BREAST CANCER CELL LINES**

Seyma Turkseven⁽¹⁾ - **Nikolaos Koundouros**⁽²⁾ - **Renata Soares**⁽¹⁾ - **Luisa Doria**⁽¹⁾ - **Alvaro Perdones-Montero**⁽¹⁾ - **Simon J. Cameron**⁽¹⁾ - **George Pouligiannis**⁽²⁾ - **Zoltán Takáts**⁽¹⁾

⁽¹⁾ University, Imperial College London / Department of Surgery and Cancer, London - ⁽²⁾ Institute, Institute of Cancer Research / Division of Cancer Biology, London

Keywords: [U-¹³C]-palmitic acid, LA-REIMS, metabolic flux analysis, breast cancer, lipogenesis

Introduction

Understanding alterations in the lipid metabolism of cancer cells is important due to the role that lipids play in membrane fluidity, energy storage and signalling processes. Using U-¹³C-labelled palmitic acid, biosynthetically incorporated into fatty acids, provides a further understanding of lipogenesis in *in vitro* studies. Here, laser assisted rapid evaporative ionisation mass spectrometry (LA-REIMS) was performed to monitor metabolic flux during cell proliferation.

Methods

SkBr3 cells were grown in DMEM containing 1% fatty acid free-BSA, and treated with [U-¹³C]-

palmitic acid under normoxic conditions for 24h. After centrifugation, cell biomass pellets were analysed directly with LA-REIMS analysis without further sample processing steps. LA-REIMS spectra were obtained in the mass range m/z 50-2500 in positive and negative ion modes and putative identification was performed using the LIPID MAPS database

Results

LA-REIMS analysis of cell lines resulted in specific lipidomic profiles that enabled the identification of different cell lines. Common adduct formations such as $[M-H]^-$ and $[M-H_2O]^-$ in negative mode and $[M+H-H_2O]^+$, $[M+H]^+$ or $[M+Na]^+$ in positive mode were detected in the metabolic flux experiment for SkBr3 cells. $[U-^{13}C]$ -palmitic acid was detected with m/z 271.2 in negative ion mode in SkBr3 cells. From the data obtained, the labelled lipid composition of SkBr3 cells in negative mode was observed to consist of phosphatidylethanolamines (e.g. PE (P-38:6) at m/z 746.5), sphingomyelins (e.g. SM (34:1) at m/z 701.5), and phosphatidylserines (e.g. PS(36:0) at m/z 790.5). In positive mode, diacylglycerides (e.g. DAG(32:1) at m/z 549.4) and ceramides (e.g. Cer(42:0) at m/z 634.6) were found in the cell line material

Conclusions

The use of stable-isotopically labelled palmitic acid in conjunction with LA-REIMS analysis allowed specific biochemical information to be collected. This approach will help to further understand how palmitic acid is being incorporated into the lipid biosynthesis pathway in a rapid and sample preparation free approach. This preliminary data demonstrates the applicability of LA-REIMS for the characterization of metabolic flux experiments using labelled metabolites.

Novel Aspect

LA-REIMS analysis can be used in metabolic flux analysis experiments in rapidly proliferating cells to model breast cancer lipogenesis – potentially identifying novel diagnostic markers and drug targets.

ThOr.06-5 RELATIVE QUANTIFICATION OF PHOSPHOLIPID SN-ISOMERS USING POSITIVE DOUBLY CHARGED LIPID-METAL ION COMPLEXES

Sven Heiles - Simon Becher - Patrick Esch

Justus Liebig University Giessen/ Analytical Chemistry, Giessen

Keywords: UVPD, FTMS, Lipidomics, Tandem MS

Introduction

Phospholipids are ubiquitous in living organisms and are involved in a plethora of biochemical processes.^[1] For glycerophospholipids assignment of double bond positions and *sn* fatty acid isomers is often not possible with routine tandem MS methods. Especially untargeted assignment of *sn* isomers is challenging. Here, we present a simple tandem MS strategy that allows phospholipid *sn* isomer assignment by phospholipid metal ion complex fragmentation.

Methods

All experiments were performed using a LTQ-FT Ultra or a Q Exactive (Thermo Fisher Scientific) mass spectrometer equipped with heated ESI sources. An in-house developed 213 nm UVPD setup^[2], CID and HCD were used for tandem MS experiments. Lipid extracts were prepared according to Schwudke et al.^[3] and phospholipase A₂ digestion was performed as described by Ekroos et al.^[4] Samples were mixed in alcohol/water and corresponding complexes were produced by ESI.

Results

Abundant lipid-metal ion complexes (LM^{2+}) complex formation was observed for phosphatidylcholines (PCs) and sphingomyelins (SMs) but not for other lipid classes. Upon LM^{2+} activation by CID, HCD or UVPD, lipid head group and fatty acid ester bond dissociation were detected. The relative abundances of the fragment ions depended on the metal ion identity. Only for iron(II) LM^{2+} complexes, head group dissociation was suppressed and only fatty acid ester bond dissociation occurred. One singly charged ion containing the dissociated deprotonated fatty acid together with Fe^{2+} ($[FA-H+Fe]^+$) and the complementary singly charged PC fragment ion were observed. For PCs containing two different fatty acid moieties, activation of LM^{2+} resulted in preferential dissociation of FAs located at *sn*2-position. PLA_2 digestion revealed that fragment ions upon CID and UVPD are *sn*2-specific but for HCD *sn*-scrambling occurs. LM^{2+} ions were utilized in MS^2 CID or UVPD experiments to quantify the relative abundance of PC *sn*-isomers in mouse organ lipid extracts.

Conclusions

Abundant ion adduction of Fe(II) to neutral PCs in ESI mass spectrometry was observed and optimized. Tandem MS experiments of these LM^{2+} complexes allow to assign fatty acid moieties and to assign *sn*-isomer in one-step tandem MS experiments. Under low energy activation conditions, relative quantification of *sn*-isomers is possible and this method can readily be applied to complex lipid mixture in a shotgun approach and potentially be coupled to LC-MS experiments.

Novel Aspect

Simple MS^2 method to quantify the relative abundance of phosphatidylcholine *sn*-isomers in complex lipid mixtures using positive ion mode.

References

1. Han X., Nat. Rev. Endocrinol., 12, 668 (2016).
2. Becher S., Spengler B., Heiles S., Eur. J. Mass Spectrom., 24 54 (2018).
3. Matyash V., Liebisch G., Kurzchalia T. V., Shevchenko A., Schwudke D., J. Lipid Res., 49 1137 (2008).
4. Ekroos K., Ejsing C. S., Bahr U., Karas M., Simons K., Shevchenko A., J. Lipid Res., 44, 2181 (2003).

ThOr.07 - Developments and applications in ion mobility MS

Chairs: Perdita Barran, Leopoldo Ceraulo

ThOr.07-1 **Keynote:** IMS-MS AS A MEANS OF REVEALING NEW STATES DURING THE MELTING OF PROTEINS

David Clemmer - Tarick El Baba - Daniel Woodall - Shannon Raab

Department of Chemistry, Indiana University

Introduction

One of the most challenging problems in biochemistry involves understanding the structures of non-native proteins that are sometimes intermediates during protein folding and unfolding events. After more than 50 years of work, experimental characterization of protein folding and unfolding usually leads to results which are described as a cooperative, two phase, transition between the folded and unfolded states – i.e., the protein appears to melt and no information about non-native states is obtained

Methods

Here we present new data from an IMS-MS analysis of simple proteins that are electrosprayed from a temperature controlled source. In this approach the source temperature is varied and IMS-MS is used to read out changes in conformations as a function of temperature. Droplets undergo an evaporative cooling phenomenon during shrinkage and this appears to trap solution structures that are not easily observed with other techniques. The IMS distributions change as the shapes of protein conformers change as the MS charge state distributions.

Results

The results suggest that the cooperative two state behavior involves other states that are captured in the IMS-MS analysis. In some examples, we find evidence for at least 10 structures that arise at slightly different transition temperatures. To date we have investigated several peptides, proteins, and protein complexes.

Conclusions

The ability to experimentally capture information about new states that are involved in folding and unfolding events may help guide theoretical efforts to model folding and unfolding processes.

Novel Aspect

IMS-MS techniques are used to capture information about non-native protein populations and structures from solution.

ThOr.07-2 **PROBING PEPTIDE CONFORMATIONAL REARRANGEMENTS WITH HIGH RESOLUTION SLIM TRAVELING WAVE ION MOBILITY-MS**

Roza Wojcik⁽¹⁾ - ***Ian K Webb***⁽¹⁾ - ***Yehia M. Ibrahim***⁽¹⁾ - ***Sandilya B.V. Garimella***⁽¹⁾ - ***Ryan Renslow***⁽²⁾ - ***Sean Colby***⁽³⁾ - ***Dennis G. Thomas***⁽³⁾ - ***Richard D. Smith***⁽¹⁾

⁽¹⁾ PNNL, Integrative Omics, Richland - ⁽²⁾ PNNL, Chem Bio & Exposure Sci Team, Richland - ⁽³⁾ PNNL, Computational Biology, Richland

Keywords: High resolution ion mobility, SLIM, peptides

Introduction

Multiply charged peptide ions in the gas phase can adopt distinct conformations, trapped in local energy minima [1]. The relative stability of these conformations provides insight into the types of secondary structures adopted by peptide ions and the mechanisms of their conformational rearrangements, such as proline cis/trans isomerization and intramolecular proton transfer. High resolution ion mobility (IM) separations are useful for probing peptide gas phase secondary structures, as they can potentially resolve multiple kinetically trapped conformers of peptides with low activation barriers for rearrangement. In this study we have employed high resolution SLIM traveling wave (TW) IM separations to study kinetically trapped peptide conformations. We investigated the impact of various IM parameters on peptide structural rearrangements during storage and separations in a SLIM IM-MS platform.

Methods

Preliminary experiments were performed using 1 meter and 13 meter Serpentine Ultra-long Path with Extended Routing (SUPER) SLIM platforms [2,3]. The SLIM utilized 6-5 (6 RF strips interspersed with 5 TW segmented arrays) or similar 4-3 electrode arrangements. Samples were introduced via a nano-electrospray ion source and trapped in the funnel trap (IFT) prior to introduction to the SLIM IM module. Afterwards ions were refocused in the rear ion funnel and passed through a segmented quadrupole to the Agilent TOF instrument. IM separations in

Nitrogen and Helium were performed in the 2.5 to 4 Torr pressure range. Simion simulations of ion trajectories through a SLIM IM module were performed with Statistical Diffusion Simulation (SDS) Model.

Results

We tracked the evolution of the kinetically trapped conformations of Substance P (M+3H)³⁺ [1] and other peptide standards as a function of IM parameters, such as the TW electric field as well as source conditions, such as RF ion funnel amplitudes. We estimated the extent of collisional heating in the SLIM platform by monitoring ion losses and fragmentation during storage and separation of selected peptide standards. Velocity distributions of peptides during SLIM IM separations have been modeled with Simion simulations. In preliminary experiments we successfully detected multiple kinetically trapped conformations of multiply charged peptide standards, including Substance P (M+3H)³⁺, Bradykinin (M+3H)³⁺ and Melittin (M+5H)⁵⁺, using both 1 meter and longer path SLIM separations. Impact of average electric field and guard bias voltage on modelled velocity distributions of Substance P (M+3H)³⁺ correlates with the distribution of the peptide conformations in local and global energy minima in the IM experiments. Collisional cross sections of the peptides under a select range of IM parameters are in close agreement with the values experimentally determined in drift tube ion mobility [4]. We detected no ion losses for Leucine Enkephalin with over 100 meter SLIM separations as well as during storage in the SLIM platform. In future experiments we plan to use different peptide ions as well as peptide noncovalent complexes to investigate separations and possible ion activation with the use of on SLIM trapping and in other novel SLIM modules. We will compare IM experiments with molecular modeling to probe the mechanisms of peptide rearrangement during IM separations.

Conclusions

SLIM TW IM enables high resolution separations of kinetically trapped conformations whose cross sections are in close agreement with the values determined by drift tube ion mobility. It also enables separation and storage of peptides with minimized dissociation. Collisional heating the SLIM platform can be controlled using certain SLIM IM parameters, enabling studies of peptide structural dynamics and the relative stability of secondary structures during IM separations.

Novel Aspect

High resolution SLIM TW IM separations of kinetically trapped peptide conformations and the characterization of collisional heating during storage and separation in SLIM.

References

1. Silveira, J.A., et al., *From Solution to the Gas Phase: Stepwise Dehydration and Kinetic Trapping of Substance P Reveals the Origin of Peptide Conformations*. Journal of the American Chemical Society, **135**,. 19147-19153 (2013)
2. Hamid, A.M., et al., *Achieving High Resolution Ion Mobility Separations Using Traveling Waves in Compact Multiturn Structures for Lossless Ion Manipulations*. Analytical Chemistry, **88**, 8949-8956 (2016)
3. Deng, L.L., et al., *Ultra-High Resolution Ion Mobility Separations Utilizing Traveling Waves in a 13 m Serpentine Path Length Structures for Lossless Ion Manipulations Module*. Analytical Chemistry, **88**, 8957-8964 (2016)
4. Stow, S.M., et al., *An Interlaboratory Evaluation of Drift Tube Ion Mobility-Mass Spectrometry Collision Cross Section Measurements*. Analytical Chemistry, **89**, 9048-9055 (2017)

ThOr.07-3 ION MOBILITY-DERIVED COLLISION CROSS SECTION FOR MYCOTOXINS: AN INTERLABORATORY STUDY

Laura Righetti⁽¹⁾ - **Lee Gethings**⁽²⁾ - **Sara Stead**⁽²⁾ - **Chiara Dall'Asta**⁽¹⁾

⁽¹⁾ Department of Food and Drug, University of Parma, Parco Area delle Scienze 95/A - 43124 Parma (Italy) -

⁽²⁾ Waters Corporation, Wilmslow, Wilmslow

Keywords: Travelling wave ion mobility spectrometry; CCS database; Interlaboratory comparison; Mycotoxins.

Introduction

Ion mobility mass spectrometry (IM-MS) is spreading rapidly as a promising technique for simultaneous screening, quantification and identification of a large number of food contaminants. Collision cross section (CCS) measurements resulting from ion IM-MS experiments provide a promising orthogonal dimension of structural information in MS based analytical separations. As with any molecular identifier, interlaboratory standardization must precede broad range integration into analytical workflows [1].

Methods

More than 60 mycotoxins were analyzed both in standard solution and in matrix to demonstrate the reproducibility of this new molecular identifier. The reproducibility of CCS values was evaluated across two traveling wave ion mobility mass spectrometers (Synapt G2 HDMS and Vion IMS QTOF) located in different laboratories.

Results

A database containing TWIMS-derived CCS values of parent and modified mycotoxins is hereby presented. CCS measurements showed high intra- and inter-day repeatability (RSD <1.8%), they were not affected by the complexity of the investigated matrices (i.e. malt, corn flakes, maize feed, wheat) ($\Delta\text{CCS} \leq 2\%$) and they resulted independent from the concentration [2]. In addition, they were not influenced by instrument settings. CCS values derived from Synapt G2 HDMS and Vion IMS QTOF mass spectrometers were highly reproducible with interlaboratory RSD <5%. Furthermore, the novel geometry of the Vion system allows to elucidate fragmentations that might occur for labile mycotoxins during different stages of the ESI-TWIMS QTOF analysis. Indeed, it is possible to state if dissociation occurs before (fragments would present shorter drift times than precursor ion) or after the TWIMS cell (drift time of the fragment ions would match that of the precursor ion). This is useful to understand which parameter has to be properly set to avoid/promote analyte fragmentation.

Conclusions

Overall, our study demonstrates the reproducibility of CCS across different instrumentations and analytical conditions. The addition of CCS in the mycotoxins screening workflows improves the accuracy and precision of analysis and would provide a higher degree of confidence in the identification workflow compared to traditional analytical approaches. Thus, many false-positive detections could be avoided by applying a maximum ΔCCS tolerance of 2%.

Novel Aspect

The low RSD observed in this interlaboratory study illustrate the potential of CCS as a molecular identifier that can be used alongside the traditional molecular identifiers.

References

1. Paglia G., Angel P., Williams J.P., Richardson K., Olivos H.J., Thompson J.W., Menikarachchi L., Lai S., Walsh C., Moseley A., Plumb R.S., Grant D.F., Palsson B.O., Langridge J., Geromanos S., Astarita G., *Analytical Chemistry*, 87, 1137-1144(2015).
2. Righetti L., Bergmann A., Galaverna G., Rolfsson O., Paglia G., Dall'Asta C., *Analytica Chimica Acta*, 1014, 50-57 (2018).

ThOr.07-4 PRACTICAL TWO-DIMENSIONAL HIGH-RESOLUTION ION MOBILITY SEPARATIONS USING FAIMS INTEGRATED WITH TRAPPED IMS

Kevin Jeanne Dit Fouque⁽¹⁾ - **Jacob Porter**⁽¹⁾ - **Francisco Fernandez-Lima**⁽¹⁾ - **Alexandre Shvartsburg**⁽²⁾

⁽¹⁾ Florida International University, Chemistry, Miami - ⁽²⁾ Wichita State University, Chemistry, Wichita

Keywords: Ion Mobility, Differential IMS, Trapped IMS, Time-of-Flight MS, Histones

Introduction

Ion Mobility Spectrometry (IMS) has emerged as a powerful broadly useful approach in conjunction with mass spectrometry, replacing or complementing condensed-phase separations. However, the resolving power of IMS remains short of that delivered by best liquid chromatography (LC) methods. The key step to frontline LC today has been transitioning from 1-D to 2-D separations that combine significantly orthogonal LC techniques for superior overall peak capacity.

Methods

While linear IMS relies on the absolute mobility, field asymmetric waveform IMS (FAIMS) is based on its increment in strong fields. As the values of a function and its derivative are *a priori* independent, these dimensions can be orthogonal and FAIMS/IMS resolves more features than either stage alone [1 - 3]. However, the first implementation involving drift-tube IMS required impractical FAIMS system floating [1]. We now coupled FAIMS to trapped IMS on a ToF platform.

Results

The FAIMS/TIMS system comprises a curved-gap FAIMS device that employs a bisinusoidal waveform with flexible profile and amplitude up to 5 kV. The unit was attached to the front of prototype Bruker TIMS instrument via an external holder, without breaking vacuum. Low MS inlet voltage permits operation without extremely complex electric floating of FAIMS system that has precluded its broad use with drift-tube IMS that requires high inlet potentials. The new platform was initially evaluated for tryptic digests, complete histone tails post-translationally modified on variant sites, and proteins. Substantial orthogonality between FAIMS and TIMS inferred from 1-D separations [3] was directly confirmed in the 2-D arrangement. In particular, some histone proteoforms co-eluting in FAIMS and TIMS individually were fully distinguished in the 2-D case. Of most interest are the findings for large proteins (e.g., albumin), where the electric dipole alignment in FAIMS (but not TIMS) results in strong orthogonality between the dimensions.

Conclusions

Ready integration of FAIMS with TIMS opens the door to practical 2-D FAIMS/linear IMS separations benefiting from exceptional resolving power (up to 400) available from TIMS. Effective ion focusing in curved-gap FAIMS filters provides high ion transmission, with overall signal similar to that in original TIMS/MS platform. Substantial orthogonality between the dimensions enables novel demanding applications.

Novel Aspect

Coupling of FAIMS with Trapped IMS produces 2-D ion mobility separations with peak capacity over 1,000, readily implementable on a commercial instrument.

References

1. Tang, K.; Li, F.; Shvartsburg, A.A.; Strittmatter, E. F.; Smith, R. D. *Anal. Chem.* 77, 6381 (2005).
2. Shvartsburg, A.A.; Li, F.; Tang, K.; Smith, R.D. *Anal. Chem.* 78, 3304 (2006).
3. Garabedian, A.; Baird, M.A.; Porter, J.; Jeanne Dit Fouque, K.; Shliaha, P. V.; Jensen, O. N.; Williams, T.D.; Fernandez-Lima, F.; Shvartsburg, A.A. *Anal. Chem.* 80, 2918 (2018).

ThOr.07-5 COLLISION CROSS SECTIONS OF PHOSPHORIC ACID CLUSTER ANIONS AND THEIR USE AS CALIBRANTS

Valentina Calabrese⁽¹⁾ - **Hélène Lavanant**⁽¹⁾ - **Frédéric Rosu**⁽²⁾ - **Valérie Gabelica**⁽³⁾ - **Carlos Afonso**⁽¹⁾

⁽¹⁾ Normandie Univ, UNIROUEN, INSA Rouen, CNRS, COBRA, Rouen - ⁽²⁾ Institut Européen de Chimie et Biologie (IECB), CNRS, UMS 3033, Pessac - ⁽³⁾ Univ. Bordeaux, INSERM, U869, ARNA Laboratory, Pessac

Keywords: IMS-MS, CCS, clusters, phosphoric acid

Introduction

The coupling of ion mobility and mass spectrometry (IMS-MS) is spreading as a supplemental method of structural characterization to mass spectrometry. Several IMS-MS techniques necessitate the use of reference ions and calibration to measure the collision cross section (CCS) of the but the pool of reference ions is limited. Since a first study in 2014 [1], we study phosphoric acid cluster ions in view of investigating their potential as calibrants.

Methods

Phosphoric acid solutions (5 to 50 mM) in water and acetonitrile were analyzed using electrospray ionization in the negative ion mode. Collision cross section were measured in helium drift gas using an Agilent 6560 drift tube ion mobility (DTIMS) Q-TOF instrument, equipped with an electrospray ionization source. Ion mobility mass spectrometry measurements were also carried out on a SYNAPT G2 HDMS (Waters), that operates a travelling wave ion mobility device.

Results

Collision cross sections from drift tube ion mobility spectrometry in helium gas were obtained for singly, doubly and triply charged phosphoric acid cluster anion. We studied the relation between collision cross sections and the aggregation number. Phosphoric acid clusters were modelled by non-overlapping spheres that adopt a tight three dimensional packing and collision cross sections were calculated using the projection approximation method. Using such minimal clusters of spheres as rough models for the cluster of phosphoric acid molecules, we could reproduce the growth of the CCS of the clusters by adapting the radius of the spheres that represent the phosphoric acid molecules.

Collision cross sections in helium were correlated to the collision cross sections in nitrogen, and their spherical models compared.

Both collision cross sections in helium and nitrogen were used to calibrate a travelling wave ion mobility instrument, and the resulting accuracy were compared using published reference compounds.

Conclusions

As phosphoric acid cluster anions are isotropic clusters that can be regarded as nanodroplets of a very accessible and simple compound, they could be used as universal calibrants for collision cross section measurement for ion mobility mass spectrometry.

Novel Aspect

Universal calibrants for ion mobility mass spectrometry using isotropic clusters.

References

1. Lavanant H., Tognetti V., Afonso C., J. Am. Soc. Mass. Spectrom., 25, 572 (2014).

ThOr.08 - Petroleomics, hydrocarbons and biofuels

Chairs: Alan G. Marshall, Giorgio Mellerio

ThOr.08-1 **Keynote:** THE COMPOSITIONAL AND STRUCTURAL CONTINUUM OF PETROLEUM REVEALED BY FT-ICR MASS SPECTROMETRY

Ryan P. Rodgers⁽¹⁾ - **Martha L. Chacón-Patiño**⁽¹⁾ - **Steven M. Rowland**⁽¹⁾ - **Sydney F. Niles**⁽²⁾ - **Donald F. Smith**⁽¹⁾ - **Greg T. Blakney**⁽¹⁾ - **Yuri E. Corilo**⁽¹⁾ - **Amy M. McKenna**⁽¹⁾ - **Alan G. Marshall**⁽¹⁾ - **Christopher L. Hendrickson**⁽¹⁾

⁽¹⁾ National High Magnetic Field Laboratory, Florida State University, Tallahassee, FL - ⁽²⁾ Department of Chemistry and Biochemistry, Florida State University, Tallahassee, FL

Keywords: FT-ICR, Petroleum, Petroleomics, Asphaltene, Complex Mixture

Introduction

Over the past 20 years, advances in modern ultra-high resolution mass spectrometry have forever changed the expectations of petroleum analyses. Facilitated by advances in ionization methods (that enable molecular-level analysis irrespective of initial or final boiling point), ultra-high resolution Fourier Transform Ion Cyclotron Resonance mass spectrometry (FT-ICR MS) routinely resolves and identifies tens-of-thousands of species in individual petroleum samples.

Methods

Petroleum samples were analyzed by FT-ICR mass spectrometry (operating at 9.4 and 21T) by atmospheric pressure photo-ionization (APPI) and/or (+/-) electrospray ionization (ESI) at concentrations between 100 -300 µg/mL. Data was acquired and processed with Predator and Petro-Org software platforms.

Results

High field FT-ICR MS has been used to revisit, and ultimately support a decades old theory of the composition of petroleum (the Boduszynski Continuum) and expose the compositional constraints (in class, type, and carbon number) imposed by boiling point. Such compositional information facilitates the accurate boiling point prediction of individual elemental compositions based on mass measurement alone, and is now poised to push forward to the molecular-level modeling of heavy ends (Petroleomics). The structural continuum is completed through the infrared multi-photon dissociation of mass isolated segments, and reveals both island (single aromatic core) and archipelago (multi-aromatic core) molecules. Thus, the current results strongly suggest that the decades-old Boduszynski model, which was largely ignored, is in fact correct, and the Yen-Mullins structural model, which is now largely accepted, is in fact, incorrect. Work supported by NSF Division of Materials Research through DMR-11-57490, State of Florida, Florida State University, and the Future Fuels Institute.

Conclusions

The current clarity on the compositional and structural continua of petroleum is heavily dependent on chromatography. Thus, current and future prospects in the on-line and off-line analysis of petroleum samples will be discussed. Work supported by NSF Division of Materials Research through DMR-11-57490, State of Florida, Florida State University, and the Future Fuels Institute.

Novel Aspect

Ultra-high resolution FT-ICR MS analysis of petroleum has revealed a compositional and structural continuum of species that changes predictably as a function of boiling point.

ThOr.08-2 **EVIDENCE OF ISLAND AND ARCHIPELAGO STRUCTURES IN DISTILLABLE AND TRULY NON-DISTILLABLE FRACTIONS OF VACUUM RESIDUES**

Diana Catalina Palacio Lozano⁽¹⁾ - ***Alexander Guzman Monsalve***⁽²⁾ - ***Enrique Mejía Ospino***⁽³⁾ - ***Mark P. Barrow***⁽¹⁾

⁽¹⁾ University of Warwick, Chemistry, Coventry - ⁽²⁾ Ecopetrol, Instituto Colombiano del Petróleo, Piedecuesta - ⁽³⁾ Universidad Industrial de Santander, Chemistry, Bucaramanga

Keywords: island, asphaltenes, archipelago, non-distillable, vacuum residues, dissociation

Introduction

Vacuum residues (VR) can be separated by molecular distillation into distillable fractions and a truly non-distillable fraction with boiling points above 687 °C. In this work, infrared multiphoton dissociation (IRMPD) was coupled to Fourier transform ion cyclotron mass spectrometry (FTICR MS) and used to determine the aromatic cores of compounds found within the vacuum residues, as a function of boiling point.

Methods

Molecular distillate fractions from two VRs were analysed, with boiling points of: 540-603 °C, 603-645 °C, 645-687 °C and, finally, 687+°C (non-distillable), all in atmospheric equivalent temperature (AET). MS/MS spectra were acquired using atmospheric pressure photoionization (APPI) coupled to a 12 T Bruker solariX FTICR MS (Bruker Daltonik GmbH, Bremen, Germany) and IRMPD (25 W CO₂ laser).

Results

Precursor ions were externally accumulated for a period of 0.5- 2 s prior to being transferred to the ion cyclotron resonance cell. It was found that the number of peaks per dalton detected increases exponentially as the boiling point of the fraction increases. A depletion of the relative intensity of polycyclic aromatic hydrocarbons and low heteroatomic-containing species was also observed as the boiling point increased, as expected using Boduszynski's model. [1] Species with up to four heteroatoms can be found in the non-distillable sample. The class distribution of the fragments indicates that up to two heteroatoms can be part of the aromatic core of the VR species, indicating that some heteroatoms can be present in the alkyl chains. MS/MS results also indicate that structural motifs comprise mainly island type structures for the distillable samples. In contrast, structures with two or more fused aromatic rings connected with alkyl chains (archipelago) are the main constituents of the non-distillable fractions of VRs.

Conclusions

Evidence of island and archipelago motifs were found in VR constituents. According with the number of double bond equivalents (DBE) of the fragments, structures of a vacuum residue can have aromatic cores with up to 8 fused benzene rings and up to two heteroatoms. Archipelago

structures are concentrated at the heavy ends fractions of crude oils, in particular in the truly non-distillable fractions.

Novel Aspect

Use of MS/MS to predict the size of the aromatic core and observation that distillable fractions contain mainly island structures.

References

1. K. H. Altgelt and M. M. Boduszynski. Taylor & Francis, 1993.

ThOr.08-3 DIRECT INLET PROBE - HIGH RESOLUTION TIME-OF-FLIGHT MASS SPECTROMETRY FOR THE DESCRIPTION OF HIGH BOILING PETROLEUM FRACTIONS

Uwe Käfer⁽¹⁾ - ***Mohammad Saraji***⁽²⁾ - ***Jürgen Wendt***⁽³⁾ - ***Thomas Wilharm***⁽⁴⁾ - ***Ralf Zimmermann***⁽⁵⁾ - ***Thomas Gröger***⁽¹⁾

⁽¹⁾ Helmholtz Center Munich, CMA, Munich - ⁽²⁾ Photonion, R&D, Schwerin - ⁽³⁾ LECO Europe, European application center, Berlin - ⁽⁴⁾ ASG, Analytik Service Gesellschaft, Neusäss - ⁽⁵⁾ University of Rostock, Department for analytical chemistry, Rostock

Keywords: Bitumen, crude oil, direct inlet probe, high resolution mass spectrometry, ionization

Introduction

High boiling petroleum fractions remain challenging to analyze on a molecular level due to their complexity and physical properties, which make them unsuitable for chromatographic techniques. Therefore, High resolution mass spectrometry has become a key technique for the analysis of these complex materials.

In this study, we used a direct inlet probe (DIP) coupled to High resolution time of flight mass spectrometry for detailed chemical characterization.

Methods

Crude oil vacuum residues and bitumen were analyzed with a multi-reflection-tof mass spectrometer (LECO Pegasus HRT) equipped with a direct inlet probe. To volatilize the sample components, a temperature gradient from 50°C to 400°C was applied. Beside electron impact ionization (EI), two soft ionization techniques, namely single photon ionization (SPI) and chemical ionization (CI) were utilized on the same system.

Results

DIP-HRMS allowed the introduction of petroleum residue samples directly into the ion source, where the temperature was linearly increased via a temperature program. Substances could be separated according to volatility and classified into chemical groups by high resolution mass spectrometry. The decreased pressure in the ion source reduces the effective boiling point, which enabled the destruction-free vaporization of high boiling substances. With the time resolved high resolution mass spectra, elemental compositions could be assigned, and homologues series of typical petro-chemical compound classes could be identified. Moreover, several bitumen samples of different origin and varying physical properties could be discriminated on a molecular level. Preliminary results showed that CI and SPI can be advantageous as complementary ionization techniques to EI because of a reduced fragmentation tendency and different selectivities.

Conclusions

Summarized, DIP-HRMS showed to be a powerful technique for the detailed description of complex, high boiling matrices. The combination of thermal separation and exact mass

information enabled identification of characteristic molecular structures. Moreover, it was possible to use the chemical fingerprint to discriminate between vacuum residues which were derived from different crude oils and refineries.

Novel Aspect

We present a new approach for in-depth analysis of high boiling samples based on the combination of thermal analysis in vacuum and high resolution mass spectrometry

ThOr.08-4 CHARACTERIZATION OF ASPHALTENES PRECIPITATED AT DIFFERENT SOLVENT POWER CONDITIONS USING MAGNETIC RESONANCE MASS SPECTROMETRY (MRMS)

Matthias Witt⁽¹⁾ - **Markus Godejohann**⁽²⁾ - **Sven Oltmanns**⁽³⁾ - **Estrella Rogel**⁽⁴⁾

⁽¹⁾ Bruker Daltonik GmbH, MRMS Solutions, Bremen - ⁽²⁾ Bruker BioSpin GmbH, Hyphenation AIC, Rheinstetten - ⁽³⁾ Bruker Optik GmbH, Applications, Bremen - ⁽⁴⁾ Chevron, Downstream Technology and Services Department, Downstream

Keywords: Asphaltene, Petroleomics, Complex Mixture

Introduction

Asphaltene precipitation is affected by factors such as heptane to crude oil ratio (HCOR). This means that for the same crude oil, a different set of molecules would precipitate forming a distinct asphaltene depending on these factors. It is expected that molecules more prone to precipitate are going to be those with the higher aromaticity, higher molecular weight and/or higher polar nature.¹⁻² Effect of solvent power on asphaltene properties HCOR is presented.

Methods

Asphaltenes were separated from a medium crude oil (API of 30.4). Different HCOR were used to induce precipitation. Asphaltenes were analyzed by MRMS. Samples were analyzed using APPI and LDI in positive ion mode with resolving power of 1.300.000 at m/z 400. Analysis of data including calculation of molecular formulas and relative abundances of compound classes was performed using PetroOrg 10.0 (Florida State University).

Results

Using LDI it was observed that HC classes as well as classes containing one heteroatom became less abundant as the heptane to crude oil ratio increases (decrease of solvent power). The opposite happened for most classes containing two or more heteroatoms. Similar results were found when positive APPI was used. Based on weighted relative abundances, DBE distributions for the samples were shifted to larger DBEs with higher solvent ratios (lower solvent power) using APPI. In contrast, molecular size distribution did not show a clear tendency. Calculations indicated that while differences in average DBE can be as large as 6 units between different HCOR, the maximum difference in molecular size is around 1.5 carbon atoms. These tendencies were also observed for individual classes for the class HC using LDI. The appearance of molecules with larger DBEs as the solvent power of the solvent decreases leads to an increase in the aromaticity. These results were in agreement with elemental analysis as well as IF-IR and ¹H-NMR results.

Conclusions

Composition of the asphaltenes changes dramatically depending on the solvent power of HCOR. Decreases in solvent power produce asphaltenes that are more aromatic. Asphaltenes become enriched in molecules with larger DBEs and with more heteroatoms per molecule. No significant changes in molecular size were observed as a function of HCOR. FT-IR and ¹H-NMR measurements confirm the tendencies evaluated by MRMS.

Novel Aspect

APPI and LDI MRMS can be used to study aromaticity and class distribution of precipitated asphaltenes at different solvent power.

References

1. Andersen, S. I. et al., *Fuel Sci. Technol. Int.*, 12, 51-74 (1994).
2. Carbognani Ortega, L. et al., *Energy Fuels*, 29, 3664-3674 (2015).

ThOr.08-5 DETAILED MOLECULAR CHARACTERIZATION OF THE UPGRADING PROCESS OF PYROLYSIS OIL BY ULTRAHIGH RESOLUTION MASS SPECTROMETRY

Yun Xu⁽¹⁾ - Zhengwen Cao⁽²⁾ - Ferdi Schüth⁽²⁾ - Wolfgang Schrader⁽¹⁾

⁽¹⁾ Max-Planck Institut für Kohlenforschung, Mass Spectrometry, Mülheim an der Ruhr - ⁽²⁾ Max-Planck Institut für Kohlenforschung, Heterogeneous Catalysis, Mülheim an der Ruhr

Keywords: Ultrahigh-resolution mass spectrometry, pyrolysis oil; biofuel upgrading

Introduction

Pyrolysis oil from biomass is a potential alternative substitution for fossil fuel. A significant drawback of pyrolysis oil is that high acidity, instability and oxygen content limits its usage as transport fuel. Much work has been done to improve the quality of pyrolysis oil. One effective way is hydrotreating approach which helps to remove oxygen and store hydrogen as well to produce petro-like fuel.

Methods

Organosolv lignin was used for pyrolysis with the condition of argon atmosphere, 600 °C at the rate of 80 °C min⁻¹. Detailed mass spectrometric measurements of the reactions after hydrotreating were carried out using a prototype-Orbitrap Elite equipped with an APCI source.

Results

Detailed ultrahigh resolution mass spectrometric analysis shows to be an excellent method to understand the detailed conversion reaction of biomass into a fuel. Here, the detailed molecular transformation can be observed at different steps of the reaction. These steps involve the pyrolysis of lignin into a first product, which then is catalytically transformed into a biofuel under hydrogen atmosphere. The results show that depending on the amount of hydrogen present during the hydrotreating process, different types of reaction products can be observed with varying degree of aromaticity. At the end, the hydrotreating products are composed of mainly hydrocarbon species and some minor oxygen-containing compounds, containing up to only four oxygens per molecule. The formation of lower DBE compositions coincides with disappearance of higher DBE compositions. And this trend can also be indicated by comparing the calculated median DBE/C values.

Conclusions

Catalytic hydrotreating from lignin derived pyrolysis oil can be fully observed using ultrahigh resolution mass spectrometry. Detailed observation of the different reaction steps allows better understanding of the conversion process.

Novel Aspect

The catalytic hydrotreating process of pyrolysis oil upgrading can be observed in detail to allow better understanding of chemical process.

ThOr.09-1 **Keynote: BIOIMAGING OF NANOPARTICLES IN CELLS AND TISSUE BY LASER ABLATION INDUCTIVELY COUPLED PLASMA MASS SPECTROMETRY**

Norbert Jakubowski - Heike Traub

Bundesanstalt für Materialforschung und -prüfung (BAM), Inorganic Trace Analysis, Berlin

Keywords: metallic nanoparticles, imaging, biological cells, tissues, inorganic mass spectrometry

Introduction

The interaction of nanoparticles (NPs) with cells and tissues has become a growing field of interest in consumer care and safety, nanotoxicology and medical applications e.g. as the ranosticor diagnostic agents. It will be shown that laser ablation inductively coupled plasma mass spectrometry (ICP-MS) is a new analytical tool which can visualize the distribution of NP and which provides quantitative data at cellular levels.

Methods

In our work we have used laser ablation coupled to ICP-MS in a differential line scanning mode to study and image the up-take and distribution of metallic nanoparticles (Au, Ag, TiO₂) in single cells. Different types of adherent cells have been grown on microscopic slides and were incubated with different types of metallic nanoparticles at different time frames and concentration.

Results

We have studied the effect of concentration, incubation time, nanoparticle diameter and surface coatings in fibroblast cells. After development of a calibration strategy using nanoparticle suspensions of known particle concentrations we could convert integrated intensities measured for a cell into the number of particles per cell. A similar strategy has been applied to convert measured intensities in tissue samples for Fe₂O₃-nanoparticles applied as contrast agents for magnetic resonance imaging (MRI) to enhance the signals in atherosclerotic plaques. By doping the oxide particles during synthesis by lanthanide elements we could easily differentiate and quantify natural (endogenous) from exogenous iron in the tissue.

Conclusions

We have optimized a laser ablation system coupled to ICP-MS to achieve sub-cellular resolution. By application of a matrix matched calibration method we can convert the measured intensities of NPs directly into the number of particles up taken by cells. The sensitivity is sufficient to detect medium sized NP as single particle or smaller once (<20 nm) as agglomerates in cells and tissue.

Novel Aspect

We discuss the concept and present the technical features of a novel quantitative elemental microscope based on inorganic mass spectrometry.

ThOr.09-2 **NANOSTRUCTURED TUNGSTEN OXIDE SUBSTRATE WITH OXYGEN VACANCIES FOR EFFICIENT SURFACE-ASSISTED LASER DESORPTION/IONIZATION MASS SPECTROMETRY ANALYSIS**

Qiang Ma

Chinese Academy of Inspection and Quarantine, Beijing

Keywords: Surface-Assisted Laser Desorption/Ionization Mass Spectrometry; Nanostructured Tungsten Oxide Substrate with Oxygen Vacancies

Introduction

The development of matrix-assisted laser desorption/ionization mass spectrometry (MALDI-MS) for characterizing molecules was a significant breakthrough that build upon the 2002 Nobel-prize winning studies of Tanaka on soft desorption ionization methods for mass spectrometric analysis. While the vacuum environment and matrix preparation make the method time-consuming and laborious, ambient ionization mass spectrometry is becoming increasingly attractive because of its simple sample pretreatment process. More than 40 ambient ionization methods have been developed since 2004 and surface-assisted laser desorption/ionization mass spectrometry (SALDI-MS) is a very important branch in the ambient MS analysis. The search for more sensitive, versatile, and robust substrate for SALDI-MS becomes particularly important.

Methods

Sea urchin-like $W_{18}O_{49}$ nanoparticles was synthesized as follows: WCl_6 dissolved in ethanol followed by hydrothermal reaction at 180 °C for 24h, and a blue flocculent precipitate was collected and diluted with ethanol to obtain the suspension. All MS analysis were performed using a home-made SALDI-MS system, and capillary nanoESI was used as a post-ionization method. For each analysis, drop a spot of $W_{18}O_{49}$ suspension on the aluminum plate as substrate, and then 5 μ L of sample solutions. The laser used was a 980 nm continuous-wave laser without lens for focusing. The relative distance and angel between components undertook a systematic optimization to acquire satisfactory MS signals.

Results

Sea urchin-like $W_{18}O_{49}$ nanoparticles have been successfully prepared via hydrothermal process. A very large absorption tail present in the visible and near infrared regions of the absorption spectrum gives clear evidence that the $W_{18}O_{49}$ nanoparticles contain a large number of oxygen vacancies, which become less for $W_{17}O_{47}$ and even be missing for WO_3 . The characteristic structure showed that compared with other nanomaterials, this substrate has a maximum energy absorption from 980 nm laser and function as “energy converter”, realizing efficient desorption of analyte. For MS analysis, the “pulse” form total ion chromatogram signal profile occurred since the power of the laser came to 4 W/cm^2 . This characteristic property gives much higher MS response and better S/N than paper spray (PS) and nanoESI analysis of 5 μ L, 20 ppb diethylhexyl phthalate (DEHP) solution in positive mode and perfluorooctane sulphonate (PFOS) in negative mode, respectively. We also noticed that there was about a 60-s delay between the signal excitation and laser irradiation, and it became shorter as the laser power increased. When the power came to 7 W/cm^2 , the delay was less than 10 s. The photothermal analysis of the $W_{18}O_{49}$ substrate dropped on the aluminum plate showed that the temperature increased faster and had a higher final value as the laser power increased. We propose that there would be a positive correlation between the signal excitation and heat accumulation but it is not a decisive factor, since the highest final temperature was about 120 °C during testing, and such a low temperature was not enough for the thermal desorption. We compared WO_3 as substrate for MS analysis and the results showed that neither oxygen vacancy nor instrumental signal. Also, for 60 min irradiation at 6 W/cm^2 , the Raman spectrum showed a certain degree of oxidation for $W_{18}O_{49}$. Corresponding MS signal for $W_{18}O_{49}$ decreased.

Novel Aspect

Novel sea urchin-like $W_{18}O_{49}$ nanoparticles with oxygen vacancies was applied for SALDI-MS analysis with significantly increased desorption/ionization efficiency.

References

1. G. Xi, S. Ouyang, P. Li, J. Ye, Q. Ma, N. Su, H. Bai, C. Wang, *Angew. Chem. Int. Ed.* 124 (2012) 2445-2449.
2. J. Zhang, Z. Li, C. Zhang, B. Feng, Z. Zhou, Y. Bai, H. Liu, *Anal. Chem.* 84 (2012) 3296-3301.

ThOr.09-3 SUPRAMOLECULAR NANOREACTORS STUDIED BY ION MOBILITY MASS SPECTROMETRY

Elina Kalenius⁽¹⁾ - ***Tung-Chun Lee***⁽²⁾ - ***Perdita Barran***⁽³⁾

⁽¹⁾ University of Jyväskylä, Department of Chemistry, Jyväskylä - ⁽²⁾ University College London, Nanomaterials Chemistry Lab, London - ⁽³⁾ University of Manchester, Manchester Institute of Biotechnology, Manchester

Keywords: Supramolecular chemistry, catalysis, nanoreactors, ion mobility mass spectrometry, CID

Introduction

Cucurbiturils (CB[n]) are pumpkin shaped macrocycles of glycouril [1]. They are water soluble compounds which can host small molecules in their cavity and operate as nanoreactors [2]. Inclusion complex formation and catalytic activity of CB[n]s in dissociation reactions of azoalkanes have been studied using MS/MS and ion mobility MS experiments.

Methods

Complex formation and MS/MS have been studied by ESI-Q-TOF MS (QSTAR from ABSciex) and Agilent 6560. Azoalkane inclusion inside nanoreactor have been followed by DTIM-MS using Agilent 6560 and modified Waters Synapt. Collision cross sections (CCS) have been compared to DFT calculated structures.

Results

Complex formation between CB[n] (n=6, 7 or 8) and azoalkanes was followed by ESI-MS. The inclusion (*endo*) complex formation was verified by DTIM-MS and it mainly depended on steric factors. Instead of *endo* complexation, formation of *exo* complexes was observed with smaller CBs. Size of CB[n] (n=6, 7 or 8), size of azoalkane and protonation state of the formed complex greatly affected the products formed during thermal activation of isolated complexes. In fact, different reaction pathway resulting in different end-products was observed for azoalkane complexed in nanoreactor on +1 and +2 charge states and for uncomplexed azoalkane.

Conclusions

CB[n] based nanoreactors are not only catalyzing reactions, but change substrate reactivity and enable formation of new reaction products. DTIM-MS is a valuable tool in analysis of structural features and especially *endo*-complexation in nanoreactors.

Novel Aspect

We have recently discovered that nanoreactors don't only catalyse reactions, but also enable formation of completely new products. DTIM-MS used for structural analysis of nanoreactors.

References

1. K. I. Assaf, W. M. Nau, *Chem. Soc. Rev.*, **2015**, 44, 394–418.
2. T.-C. Lee, E. Kalenius, A. I. Lazar, K. I. Assaf, N. Kuhnert, C. H. Grün, J. Jänis, O. A. Scherman, W. M. Nau, *Nat. Chem.*, **2013**, 5, 376–382.

ThOr.09-4 H₂ DETECTION IN-SITU BY MEMBRANE INTRODUCTION MASS SPECTROMETRY

Alice Rousseau - ***Ana-Laura Luna Baron*** - ***Andrea Mendez*** - ***Joel Lemaire*** - ***Michel Heninger*** - ***Hynd Remita*** - ***Essyllt Louarn***

Keywords: MIMS, energy, hydrogen, light rare gases, real-time analysis

Introduction

Direct measurement of gases formed in hydrogen reactors based on nanomaterials [1] is often challenging as many molecules may be of interest, rapid response is necessary and small volumes are involved. We propose an instrument, based on mass spectrometry, for small inorganic gases analysis (H₂, CO₂, O₂, NO, Ar,...) specifically adapted to the hydrogen production applications.

Methods

Gas analysis in the reactor is processed by a membrane interface (MIMS) and a quadrupole mass spectrometer. As water from the permeate may interfere with the hydrogen signal [2], we added a cryotrap at -100°C before the vacuum cell. The pumping system was modified to increase compressibility of hydrogen. Calibration of H₂ signal is managed by a gas diluter. The apparatus was tested in a H₂ formation reactor.

Results

The apparatus was tested in dry and in-water conditions. Headspace analysis was not possible as signal was not stable due to droplet formation on the membrane. For in-water analysis stirring of the solution is necessary to maintain constant signal.

The modifications we made on the instrument improved the detection and stability of the H₂⁺ (*m/z* 2) signal. Signal had a linear response to increasing concentrations. Detection limit was 2 ppm. First, use of specific material for the vacuum cell decreased hydrogen and water adsorption and ultimately lowered detection limit. Second, memory effect due to low compressibility of hydrogen was solved by the use of a second turbo pump. Blanks were constant before and after introduction of high quantity of H₂ for a long period. Finally, addition of a cryotrap decreased drastically (more than 99%) water amount in the permeate, and therefore detection limit of hydrogen was greatly improved in particular in wet conditions.

Conclusions

The constructed prototype allowed for improved H₂ detection in mass spectrometer. Signal is more stable and detection is made possible directly in the water solution, with no headspace required. Concomitant analysis of O₂ and CO₂ was produced in a hydrogen formation reactor.

Novel Aspect

New instrument specifically adapted to analysis in real-time of H₂ and different other gases directly connected to a reactor to monitor hydrogen production from nanomaterials.

References

- 1 Luna-Barron A. et al., *Applied Catalysis B: Environmental*, 191, 18-28 (2016)
2. King S. and Price S., *Int. J. Mass Spectrom.*, 277,84-90 (2008)

ThOr.09-5 **DEVELOPMENT OF ISOLATION TECHNIQUES AND CAPILLARY ELECTROPHORESIS MASS SPECTROMETRY (CESI-MS) METHODS FOR THE CHARACTERISATION OF NANOMATERIAL PROTEIN CORONAS**

Klaus Faserl⁽¹⁾ - **Andrew Chetwynd**⁽²⁾ - **Iseult Lynch**⁽³⁾ - **James Thorn**⁽²⁾ - **Herbert Lindner**⁽¹⁾

⁽¹⁾ University of Innsbruck, Department of Medical Chemistry and Biochemistry, Innsbruck - ⁽²⁾ AB Sciex UK Ltd, Separations, Warrington - ⁽³⁾ University of Birmingham, School of Geography Earth and Environmental Sciences, Birmingham

Keywords: CE-MS, Protein-corona, nanoparticle, proteomics

Introduction

Since the nanoparticle corona rose to eminence a decade ago, it has been investigated by a wide range of researchers trying to characterise and derive its significance [1]. To date, LC-MS is the instrument of choice for the analysis of the protein corona [1,2]. However, CESI-MS offers an exciting prospect to increase sample throughput and separate and detect highly polar and very large peptides which may be missed by conventional nLC-MS platforms [3].

Methods

Silica nanoparticles were incubated for 1 hr at 37 °C in plasma. The on-particle digest particles were washed followed immediately by digest. In the case of the off particle digest, NP-corona complexes were boiled in SDS buffer and run on an SDS-PAGE gel prior to in-gel digestion. A Sciex CESI 8000 equipped with a 90 cm neutral capillary with 30 kV separation voltage and 2 PSI pressure was used for peptide separation, detection was with a Thermo QExactive HF mass spectrometer.

Results

This investigation began by developing a sample preparation method to isolate the corona prior to CESI-MS analysis. Here the off-particle digest method enabled the detection of 1581 peptides which correlated to 134 protein groups. Initial on-particle digest methods were limited to 1270 peptides and 147 proteins due to peptides adsorbing to the NP surface to reform a corona. To ameliorate this, the addition of urea and/or rapigest to the digest solution were investigated as potential methods to keep the digested peptides solubilised. The addition of 1.5 M urea with and without 0.1% rapigest improved both peptide and protein coverage. However, the addition of 0.1% Rapigest only enabled the detection of 2215 peptides and >200 protein groups. Furthermore, this method resulted in the least number of missed cleavages whereas any addition of urea increased the number of missed cleaves compared to the initial ammonium bicarbonate digest buffer solution.

Conclusions

The finalized sample preparation and CESI-MS methods enabled a much higher throughput analysis than a conventional nLC-MS method, 50 minutes versus 140 minutes injection to injection. This was achieved by using an on-particle digest (<12 hrs) instead of an off-particle in gel digest (1.5 days). In addition, this study introduced CESI-MS to the field of nanoparticle protein corona analysis.

Novel Aspect

This is the first example of CESI-MS being applied to the nanoparticle protein corona and the first to comprehensively compare on and off particle corona isolation methods.

References

1. Ke P., Lin S., Parak W., Davis T., Caruso F., *ACEnano*, 11, 11773-11776 (2018).
2. Carrillo-Carrion C., Carril M., Para W., *Current Opinion in Biotechnology*, 46, 106-113 (2017).
3. Chetwynd A., Guggenheim E., Briffa S., Thorn J., Lynch I., Valsami-Jones E., *Nanomaterials*, 8, 99 (2018).

ThOr.10 - Novel food; nanoparticles in food and beverage; future applications

Chair: Stefan Weigel

ThOr.10-1 **Keynote: MASS SPECTROMETRIC DETECTION AND CHARACTERISATION OF NANOPARTICLES**

Ruud Peters - Anna Undas - Martijn van der Lee

RIKILT Wageningen University & Research, Wageningen, The Netherlands

Keywords: nanoparticles, HDC, AF4, ICP-MS, spICP-MS

Introduction

Nanotechnology is on its way to become a key enabling technology and is expected to stimulate industrial growth, innovation and development in the most diverse fields. Although potential beneficial effects of nanotechnologies are generally well described, still existing uncertainties concerning environment, health and safety aspects need to be addressed to explore the full potential of this new technology. One challenge consists in the development of methods that reliably identify, characterize and quantify nanomaterials, both as substance and in various products and matrices. Mass spectrometry plays an important role here.

Methods

Recent applications of mass spectrometry for the identification and characterization of nanoparticles have been reviewed.

Results

Mass spectrometry, in the form of inductively coupled plasma mass spectrometry (ICP-MS) has been used in combination with hydrodynamic chromatography (HDC) to determine the presence and concentration of metal and metal oxide nanoparticles. This combination has been used to determine not only the presence of NP in environmental samples and food, but also to study the fate of silica nanoparticles following digestion. MALDI-TOFMS has been used to identify organic nanoparticles but can also be used to identify organic coatings on nanoparticles. In many occasions asymmetric field flow fractionation (AF4) has replaced HDC as the separation method for nanoparticles. To detect and quantify nanoparticles, AF4 is often combined with ICP-MS. AF4-ICP-MS has been used to identify and characterise nanoparticles in environmental samples and in food and personal care products. A relatively new technique is single particle ICP-MS (spICP-MS) which serves as a screening technique for nanoparticles. spICP-MS determines the particle size and the particle number- and mass-based concentration. The attractive feature is that spICP-MS determines a number-based size separation, very much in line with European recommendation for the definition of a nanoparticle. spICP-MS can also detect nanoparticles at very low concentrations which is an advantage in sample processing. Hyphenation of spICP-MS with separation techniques is also studied. The combination of HDC and AF4 with spICP-MS allows for the determination of aggregates and agglomerates since HDC and AF4 determine the hydrodynamic diameter of a particle while spICP-MS determines the spherical equivalent diameter. Very recently, time of flight mass spectrometry has been introduced in spICP-MS allowing true multi-element detection and opening new possibilities.

Conclusions

Mass spectrometry is well on its way to become an important technique in the detection and characterisation of nanoparticles in the environment and in biological and food samples. Especially spICP-MS has become an useful screening method for nanoparticles and combined with electron microscopy as a confirmation method it forms a strong combination.

Novel aspect

The identification and characterization of nanoparticles is a new area of research and one that shows many novel aspects.

ThOr.10-2 **HOW TO MANAGE MATRIX EFFECTS FOR ACCURATE NANOPARTICLE SIZING USING SINGLE-PARTICLE ICPTOFMS?**

Lyndsey Hendriks - Alexander Gundlach-Graham - Detlef Günther

ETH Zurich, Department of Chemistry and Applied Biosciences, Zurich

Keywords: nanoparticles, single-particle ICPMS (sp-ICPMS), time-of-flight MS (TOFMS), food, complex matrices

Introduction

The use of nanoparticles (NPs) in consumer products is not new. In our kitchen, inorganic NPs can be found as food additives with specific purposes: color or flavor enhancers, preservative agents and texture modifiers. Detection and characterization of NPs in these complex media is necessary to assess risk linked with their use; however, standardized protocols and methods are still lacking for NPs characterization in complex samples.

Methods

Single-Particle Inductively Coupled Plasma Mass Spectrometry (sp-ICPMS) is used to measure elemental mass of individual NPs and particle number concentrations in complex matrices. To investigate matrix effects on NP measurements, we combined ICPTOFMS (icpTOF, TOFWERK AG, Switzerland), with a dual sample introduction system, in which microdroplets doped with analytes of interest are used for calibration [1].

Results

Gold NPs were measured in diverse matrices with varying complexity, ranging from ultra-pure water to food samples (fruit juices & milk) and a buffer solution (PBS) used in biological research. In fruit juice, we observed few matrix effects—sizing of Au NPs was identical to that obtained in water. However, with Au NPs dispersed in milk or PBS, we observe an attenuation of NP and microdroplet signal intensities. Because the sensitivities of calibrant droplet and analyte NP are affected in the same way, Au NPs could be sized accurately despite this matrix effect. In contrast, without online microdroplet calibration, Au NP diameters would have been underestimated by 27% or 34% for milk and PBS matrices, respectively.

Conclusions

Microdroplets (calibrant used to determine the instrument sensitivity) are introduced concurrently with the NP-containing sample and experience the same plasma conditions. This online calibration approach can be used to compensate for plasma-related matrix effects. This setup allows for the analysis of large volume samples, rapid switching of NP-containing samples, and calibration of a wide range of NPs.

Novel Aspect

We demonstrate a successful strategy to detect and accurately size metallic NPs in various real samples with automatic compensation for plasma related matrix effects.

References

B. Ramkorun-Schmidt, S. A. Pergantis, D. Esteban-Fernandez, N. Jakubowski and D. Gunther, *Anal Chem*, 87, 8687-8694 (2015)

ThOr.10-3 **USE OF ICP-MS-BASED TECHNIQUES FOR THE CHARACTERIZATION OF NANOMATERIALS USED IN FOOD-RELATED APPLICATIONS**

Federica Aureli - Francesca Ferraris - Sara Savini - Andrea Raggi - Francesco Cubadda

Istituto Superiore di Sanità, National Institute of Health, Rome

Keywords: nanoparticles, food additives, nutrient sources, single particle inductively coupled plasma mass spectrometry, asymmetric flow field flow fractionation-UV-MALS-ICP-MS

Introduction

Applications of nanotechnologies to agricultural production, food processing, and food contact materials are rapidly developing. Such applications may bring benefits but potential risks have to be assessed and excluded since, if nanoparticles persist as such after gastrointestinal digestion, they may be absorbed in the gut. In the NANO-PERSIST project the transformations and fate of food-relevant nanoparticles during human digestion and in lysosomal conditions are being studied.

Methods

In the project, carried out in collaboration with the EC-JRC, acellular in vitro tests are applied to assess the dissolution/biopersistence of nanomaterials used in the food and nutrition sector and form the basis of the risk assessment strategy for these materials. A multi-method approach is used for the physicochemical characterization of the samples, where the results of MS-based techniques (spICP-MS, AF4-UV-MALS-ICP-MS/MS) are compared with those of centrifugal liquid sedimentation (CLS) and transmission electron microscopy (TEM).

Results

Ten different samples, belonging to four material types (namely synthetic amorphous silica, titanium dioxide, iron oxides/hydroxides, and zinc oxide) were studied. The pristine materials were characterized after dispersion with a standardized protocol. Whereas TEM was essential to study particle morphology, MS-based techniques provided excellent performances for the size characterization of the materials after proper method development. In particular for spICP-MS, the use of μ s dwell times combined with the use of ion–molecule chemistry for resolution of spectral interferences was essential to obtain quantitative data - with fast acquisition times - and with improved (lower) size LoDs.

Conclusions

ICP-MS based techniques are the most promising tools for routine characterization and detection of inorganic nanomaterials, also in complex matrices after proper sample preparation. Development of multi-technique approaches and combination of separation and sizing techniques appears to be the way forward for the physicochemical characterization of particulate materials used as food additives and nutrient sources.

Novel Aspect

For the first time a range of materials belonging to four different material types have been characterized with state-of-the-art MS-based techniques and the results compared with TEM and CLS.

ThOr.10-4 TRANSLATING INTO PRACTICE THE EFSA GUIDANCE ON RISK ASSESSMENT OF THE APPLICATION OF NANOTECHNOLOGIES IN THE FOOD CHAIN: THE ROLE OF MASS SPECTROMETRY

Federica Aureli - Francesca Ferraris - Sara Savini - Andrea Raggi - Francesco Cubadda

Istituto Superiore di Sanità, National Institute of Health, Rome

Keywords: nanomaterials, food, risk assessment, MS-based techniques, single particle inductively coupled plasma mass spectrometry

Introduction

The European Food Safety Authority has produced a Guidance on the risk assessment of nanoscience and nanotechnology applications in the food chain covering application areas such as novel foods, food contact materials and food additives. As far as analytical chemistry is concerned, the Guidance elaborates on physicochemical characterisation of nanomaterials in terms of how to establish whether a material is a nanomaterial, the key parameters that should be measured, the methods and techniques that can be used for characterisation of nanomaterials and their determination in complex matrices.

Methods

Practical examples of the key role of MS-based techniques in the implementation of this Guidance are given. In the characterisation of pristine materials, they provide information on the chemical nature of the particles to be studied, which is essential in addition to data on size and other parameters. In the determination of nanoparticles in complex matrices, owing to the possible presence of other particulate materials in the test sample, the use of a chemically specific method, e.g. relying on single particle inductively coupled plasma mass spectrometry (spICP-MS), is necessary and in many cases more viable compared to application of electron microscopy techniques, e.g. TEM-EDX.

Results

The use of spICP-MS enabled the size characterization of a pristine food-grade anatase-based titanium dioxide particulate material. Results were checked by centrifugal liquid sedimentation (CLS) and transmission electron microscopy (TEM). In the subsequent biodistribution study in mice, after oral administration of the material, ICP-MS/MS (a non-nanospecific technique) allowed to screen for the total Ti in tissues owing to the low background levels of titanium in the body. A significant deposition was ascertained in the large intestine and spICP-MS (after TMAH extraction) was then used to check the particulate nature of the Ti present. The Ti was found to be present as particles and the size distribution was determined.

Conclusions

An anatase-based titanium dioxide food additive was characterised in pristine form and its biodistribution studied after oral administration to mice according to the EFSA nano-Guidance. By using ICP-MS/MS and spICP-MS it was shown that the material is largely in the nano-range, it is absorbed in the body and it specifically accumulates in the large intestine.

Novel Aspect

The usefulness of MS-based techniques in the implementation of the novel EFSA nano-Guidance has been demonstrated.

FRIDAY, August 31st, 2018

PL06: Plenary Lecture

Chair: Simonetta Fornarini

PL06 **MASS SPECTROMETRY AND THEORETICAL CHEMISTRY IN SERVICE OF CATALYSIS RESEARCH: A MÉNAGE-À-TROIS AT ITS BEST**

Helmut Schwarz

Institut für Chemie, Technische Universität Berlin, 10623 Berlin, Germany

The ultimate goal in heterogeneous catalysis is to make use of each and every atom of supported (metal) catalysts, i.e. in the extreme to perform single-atom catalysis (SAC). While this arduous task constitutes a non-trivial, if not daunting challenge in 'real-life' chemistry, in the gas phase SAC can be achieved in a rather straightforward manner by conducting experiments with mass-selected species under (near) single-collision conditions. These mass spectrometry-based studies on isolated reactants, when complemented by state-of-the-art computational and spectroscopic work, provide an ideal arena for probing the energetics and kinetics of a chemical process in an unperturbed environment at a strictly molecular level without being obscured by ill-defined side effects. Thus, the concept of SAC can be explored or, more generally, the mechanisms of reactions and the active parts of single-site catalysts, the so-called 'aristocratic' atoms, can be identified.

Examples discussed include:

- (i) The room-temperature, cluster-oxide mediated redox reactions of the CO/N₂O couple which, arguably constitute one of the most prominent oxidation processes, and
- (ii) novel metal-mediated C–H bond activation and C–C coupling of methane, which are regarded as one of the holy grails in chemistry.

References

1. H. Schwarz, *Catal. Sci. Technol.* **2017**, 7, 4302
2. H. Schwarz, S. Shaik, J. Li, *J. Am. Chem. Soc.* **2017**, 139, 17201
3. H. Schwarz, *Angew. Chem. Int. Ed.* **2011**, 51, 1009

FOr.01 - Life sciences – Metabolomics

Chairs: Ulf Sommer, Donatella Caruso

FOr.01-1 **Keynote: METABOLOMICS AND THE EXPOSOME - APPLICATION TO CANCER EPIDEMIOLOGY**

Augustin Scalbert

International Agency for Research on Cancer (IARC), Biomarkers Group, Nutrition and Metabolism Section, Lyon, France

Keywords: *Exposome, biomarkers, metabolomics, diet, cancer risk*

The exposome, defined as all exposures to which an individual is subjected along life course, involves chemicals or metabolites that can be directly measured in biospecimens such as blood or urine. Hundreds of these chemicals or metabolites, derived from diet, drugs, pollutants and other environmental factors, can affect health and risk of diseases. High-resolution mass spectrometry (MS) techniques, combined with biostatistics and bioinformatics in metabolomic studies, are increasingly used to measure these chemicals and metabolites, in order to discover new biomarkers of exposure or new risk factors for diseases in exposome-wide association studies. These MS techniques are particularly adapted to the measurement of a large diversity of biomarkers of exposure and to assess complex exposures such as diet.

Each food is indeed characterized by a specific chemical profile, with compounds often specific for a particular food or food group that are eventually absorbed during digestion. These compounds found in the systemic circulation constitute the 'food exposome'. Examples will be presented on the identification of novel biomarkers of food intake in dietary intervention and cross-sectional epidemiological studies and on the study of their associations with cancers in recent case-control studies nested in large cohort studies, thus providing new evidence on the role of the diet in disease risk. The same metabolomic studies also provide new data on associations between exposures and the host metabolism shedding new light on mechanisms linking exposures to diseases as will be illustrated for diet, air pollution, overweight or birth

weight.

Development of these studies also faces a number of challenges. These include a more comprehensive annotation of the exposome, the measurement of components of the exposome present at low concentrations, the development of quantitative, high throughput and robust methods (targeted metabolomics) for large epidemiological studies. Recent progress in this rapidly moving field will be presented. It should contribute to the identification of novel disease risk factors and to the improvement of our understanding of disease aetiology.

FOr.01-2 TARGETED ON-LINE BREATH ANALYSIS SUPPORTS THE HYPOTHESIS OF ALTERED COLLAGEN TURNOVER IN IDIOPATHIC PULMONARY FIBROSIS

Martin Thomas Gaugg⁽¹⁾ - ***Anna Engler***⁽²⁾ - ***Pablo Martinez-Lozano Sinues***⁽³⁾ - ***Lukas Bregy***⁽¹⁾ - ***Yvonne Nussbaumer-Ochsner***⁽²⁾ - ***Lara Eiffert***⁽¹⁾ - ***Tobias Bruderer***⁽¹⁾ - ***Renato Zenobi***⁽¹⁾ - ***Malcolm Kohler***⁽²⁾

⁽¹⁾ Swiss Federal Institute of Technology, Department of Chemistry and Applied Biosciences, Zurich - ⁽²⁾ University Hospital of Zurich, Department of Pulmonology, Zurich - ⁽³⁾ University of Basel, Department of Biomedical Engineering, Basel

Keywords: on-line breath analysis, idiopathic pulmonary fibrosis, biomarker, secondary electrospray ionization – mass spectrometry

Introduction

On-line breath analysis is a powerful technique to obtain insights into the metabolism of a person. Idiopathic pulmonary fibrosis (IPF) is a chronic and poorly understood lung disease. Following a recent study that reports increased levels of collagen related amino acids in lung tissue of IPF patients using GC-MS profiling[1], we hypothesized that these altered amino acid levels might be mirrored in exhaled breath, which would allow for a non-invasive screening of IPF.

Methods

Breath analysis was performed on-line using secondary electrospray ionization-high resolution mass spectrometry (SESI-HRMS). For all previously reported target compounds where a signal was detected at the accurate mass in real-time, UHPLC-MS measurements of exhaled breath condensate (EBC) were performed, to confirm their identity. Their classification performance was estimated using 1 million leave-one-out cross-validations.

Results

On-line breath spectra were recorded for 21 IPF patients and 21 healthy controls, matched in terms of age, gender and smoking history. Within these, we could detect robust signals for proline, 4-hydroxyproline, alanine, valine, leucine/isoleucine, allysine, phenylalanine and pyroglutamic acid, most of which could also be confirmed in EBC. Six of the eight compounds showed significantly increased signal levels ($p < 0.05$) in exhaled breath of IPF patients. Additionally, we observed a strong signal correlation across subjects, and those amino acids with a higher abundance in collagen or elastin showed the most pronounced effect. This is consistent with altered collagen and elastin turnover being the underlying metabolic processes. Using the signals of all detected amino acids, we were able to obtain a cross-validated area under the receiver operating characteristic curve of 0.86.

Conclusions

Using targeted on-line breath analysis with SESI-HRMS, we could detect increased amino acid levels in IPF patients, which allowed for a good discrimination from healthy controls. This is consistent with previous metabolomic findings from lung biopsies, strongly reducing the probability of false discoveries. However, we were able to capture this information in a non-invasive and rapid fashion, underlining the strength of real-time breath analysis.

Novel Aspect

Our results suggest that on-line breath analysis using SESI-HRMS might be a valuable tool for screening of IPF, thereby reducing the time of diagnosis as well as medical costs.

References

1. Kang Y. P., Lee S. B., Lee J. M., Kim H. M., Hong J. Y., Lee W. J., Choi C. W., Shin H. K., Kim D. J., Koh E. S., Park C. S., Kwon S. W., Park S. W., *J. Prot. Res.*,15 (5), 1717-1724,2016

FOR.01-3 TARGETED METABOLOMIC ANALYSIS OF URINE FOR IMPROVING THE DIAGNOSIS OF ASTHMA AND COPD

Mona Khamis⁽¹⁾ - **Hanan Awad**⁽²⁾ - **Kevin Allen**⁽¹⁾ - **Darryl Adamko**⁽³⁾ - **Anas El-Aneed**⁽¹⁾

⁽¹⁾ University of Saskatchewan, College of Pharmacy and Nutrition, Saskatoon - ⁽²⁾ Calgary Laboratory Services, Alberta Health Services, Calgary, AB, Canada - ⁽³⁾ Department of Pediatrics, College of Medicine, University of Saskatchewan, Saskatoon, SK, Canada

Introduction

Asthma and COPD are chronic inflammatory conditions of the respiratory airways. Their differential diagnosis can be hampered by overlapping clinical presentations. Currently, physicians rely on patient history and trials of therapy. Therefore, there is a need for new diagnostic tests. A recent untargeted exploratory ¹H-NMR metabolomic study [1] suggested 50 urinary metabolites as candidate biomarkers for the differentiation of asthma and COPD.

Methods

The 50 metabolites were divided into 4 groups. Groups 1 and 2 (37 targets) contain amine and carboxylic acid-containing metabolites. Two LC-MS/MS methods for their quantification were developed using differential isotope labelling (DIL) with dansyl chloride and dimethylaminophenacyl bromide, respectively. Group 3 contains 7 metabolites of diverse structures. Therefore, HILIC-MS/MS method was developed. All methods were fully validated for clinical testing.

Results

All methods were fully validated according to the FDA and European Medicines Agency (EMA) guidelines [2, 3]. While almost all of the validation criteria were still applicable, in some situations, the guidelines were insufficient to drive the validation process. Accordingly, uncommon analytical or statistical approaches were adopted. In addition, the use of DIL with 2 mass unit difference between the derivatized analyte and its internal standards necessitates the optimization of the analytical methods to account for isotopic interferences. Non-blinded patients' samples with asthma and COPD were analyzed and the data was processed using partial least square discriminant analysis. A set of fifteen metabolites was found to differentiate asthma from COPD with high accuracy (R² 0.935 Q² 0.864). The constructed model was then tested using a blinded subset of subjects, where an accuracy of diagnosis of 80% for asthma was achieved. However, COPD had modest accuracy; examination of clinical information suggested misdiagnosis for some of these patients.

Conclusions

Novel validated analytical methods were developed for clinical application. Promising clinical results were obtained which may improve the management of airway diseases. Of the 15 metabolites identified by statistical analysis from the mass spectrometry data, eight metabolites were common, as critical for disease diagnosis, to the previously reported ¹H-NMR study [1]. Therefore, this work also validates the diagnostic importance of these metabolites. A method for group 4 metabolites is currently under development.

Novel aspect

New LC-MS targeted methods are developed, providing insights into targeted metabolomics method validation. 15 urine biomarkers are identified for respiratory illnesses

References

1. Adamko D.J., Nair P., Mayers I., Tsuyuki R.T., Regush S., Rowe B.H., *Journal of Allergy and Clinical Immunology*, 136, 571-80. e3 (2015).
2. FDA. <https://www.fda.gov/downloads/drugs/guidances/ucm368107.pdf> (2013).
3. European Medicines Agency (EMA), http://www.ema.europa.eu/docs/en_GB/document_library/Scientific_guideline/2011/08/WC500109686.pdf (2011).

FOR.01-4 **IN VIVO REAL-TIME MONITORING SYSTEM FOR METABOLITES IN A LIVING MOUSE BRAIN USING PROBE ELECTROSPRAY IONIZATION/TANDEM MASS SPECTROMETRY (PESI/MS/MS)**

Kei Zaitzu

Nagoya University, In Vivo Real-time Omics Laboratory, Institute for Advanced Research, Nagoya

Keywords: In vivo real-time monitoring, Metabolites, PESI/MS/MS, Real-time metabolomics

Introduction

Recent improvements in ambient ionization techniques have achieved real-time monitoring of biogenic compounds in organisms. Here, we developed a novel analytical system for in vivo real-time monitoring of metabolites in a living mouse brain [1]. It consists of a manipulation system and a probe electrospray ionization unit, which uses an extremely thin solid needle (tip dia.: 700 nm) for direct sampling and ionization, coupled to a tandem mass spectrometer [1-3].

Methods

The system consists of a LCMS-8040 equipped with PESI (Shimadzu) and the manipulation system (Shimadzu). To validate the system, a mouse was anesthetized and ca. 4 mm-dia. hole was drilled at the center of the skull to expose the intact brain surface. A special sample cup was fixed onto the skull, and 8 metabolites were analyzed by the system. A CB₁ agonist, which also acts as an inhibitor of TCA-cycle activity, or vehicle were administered 3 min after the analysis started.

Results

Eight metabolites including TCA-cycle intermediates were successfully monitored in real time in a living mouse brain for 3 hours. Glutamic acid levels immediately decreased after the analysis started for both CB₁ agonist-administered and vehicle mice. For CB₁ agonist-administered mouse, glucose rapidly decreased 30 min into the analysis. At the same time, citric/isocitric acid increased remarkably. ADP and TCA-cycle intermediates, succinic, fumaric and malic acids, maintained its level or weakly increased, while α -ketoglutaric acid level slightly decreased. This state continued for about 30 min, then returned to their initial levels. During these dynamics observations, glucose and citric/isocitric acids seem to fluctuate coordinately. On the other hand, these phenomena were not observed in the vehicle mouse. Additional experiments demonstrated that the extracellular glutamic acid changes monitored by microdialysis showed the opposite profile, suggesting that the present system may be capturing intracellular metabolic changes.

Conclusions

We developed a novel in vivo real-time monitoring system for a living mouse brain by combinational use of PESI/MS/MS and a manipulation unit. We applied the system to real-time

monitoring of metabolites in a living mouse brain, successfully capturing the metabolic dynamics in mouse brain over 1 hour with 20-second intervals. The system is highly expected to be expanded to real-time metabolomics in the near future.

Novel Aspect

The newly-developed analytical system achieved in vivo real-time monitoring of the metabolites in a living mouse brain.

References

1. Zaitso K, Hayashi Y, Murata T. et al., *Analytical Chemistry*, 90 (7), 4695-4701 (2018).
2. Zaitso K, Hayashi Y, Murata T. et al., *Analytical Chemistry*, 88 (7), 3556-3561 (2016).
3. Hayashi Y, Zaitso K, Murata T. et al., *Analytica Chimica Acta*, 983, 160-165 (2017).

FOr.01-5 ENERGY METABOLISM IN AGE-RELATED DISEASES: FOCUS ON MITOCHONDRIA

Nico Mitro - Donatella Caruso

Università degli Studi di Milano, Department of Pharmacological and Biomelecular Sciences, Milano

Keywords: steady-state metabolomics, fluxomics, energy metabolism, peripheral neuropathy, diabetes.

Introduction

Metabolomics is a powerful tool to gain new insights contributing to the identification of complex molecular mechanisms in human and animal cells. Our aim is to understand how metabolism is rearranged during the development of age-related diseases. In particular, we focused our attention on the role of mitochondria that represents the energy-generating hubs of the cell in the context of two pathologies related to age: the peripheral neuropathy and type 2 diabetes.

Methods

We used ESI-LC-MS/MS metabolomics of more than 500 metabolites involved in energy metabolism to investigate changes in the sciatic nerve of an experimental model of peripheral neuropathy. The same analyses were also used to characterize a novel mitochondrial regulator *in vitro* and in humans. In these latter cases, we also determined the rates of metabolic reactions by culturing cells with uniformly ¹³C-labeled nutrients (glucose, palmitate and glutamine).

Results

Mice lacking the master regulator of *de novo* lipogenesis Sterol Regulatory Element Binding Factor (SREBF-1cKO) developed peripheral neuropathy. Metabolomics indicated that SREBF-1cKO peripheral nerves experienced decreased fatty acid synthesis and glycolytic flux, but increased fatty acid catabolism and mitochondrial function. These alterations were the result of local accumulation of two endogenous PPAR α ligands, specifically identified by mass spectrometry [1].

We have identified Zinc finger CCCH-type containing 10 (Zc3h10) as novel mitochondrial regulator. Depletion of Zc3h10 results in mitochondrial dysfunction and blunted Krebs cycle flux. We have also identified a loss-of-function mutation of Zc3h10 in humans (Tyr105 to Cys105) that is associated with increased body mass index and hyperglycemia. Isolated peripheral blood mononuclear cells from Cys105 homozygotes display reduced oxygen consumption rate and decreased levels of some Krebs cycle metabolites, which all together derive in mitochondrial dysfunction [2].

Conclusions

Collectively, our results contributed to elucidate metabolic pathways modulated by two specific proteins, namely SREBF-1c and Zc3h10. Specifically, we identified changed energetic metabolites when both proteins were experimentally knock-down or functionally altered. Further, stressing the physiologic relevance of our findings, we demonstrated the key role of SREBPF-1c and Zc3h10 in peripheral neuropathy and mitochondrial dysfunction development.

Novel Aspect

Steady-state and dynamic metabolomic analyses are powerful analytical methods to unravel metabolic re-wiring in the transition from physiology to pathophysiology.

References

1. Cermenati G, Audano M, Giatti S, Carozzi V, Porretta-Serapiglia C, Pettinato E, Ferri C, D'Antonio M, De Fabiani E, Crestani M, Scurati S, Saez E, Azcoitia I, Cavaletti G, Garcia-Segura LM, Melcangi RC, Caruso D, Mitro N., *Cell Metabolism*, 21, 571-83 (2015).
2. Audano M, Pedretti S, Cermenati G, Brioschi E, Diaferia GR, Ghisletti S, Cuomo A, Bonaldi T, Salerno F, Mora M, Grigore L, Garlaschelli K, Baragetti A, Bonacina F, Catapano AL, Norata GD, Crestani M, Caruso D, Saez E, De Fabiani E, Mitro N., *EMBO Reports*, 19, pii: e45531 (2018).

FOr.02 - Informatics tools and data analysis

Chairs: Peter O'Connor, Giuseppe Avellone

FOr.02-1 **Keynote: LEVERAGING NEW ALGORITHMICS FOR BETTER MS SPECTRA**

Marc-André Delsuc

IGBMC, 1 rue laurent Fries, BP 10142, 67404 Illkirch, France

Keywords: Big-Data, Data Analysis, Isotopic Deconvolution, Pattern Matching, Compressed Sensing

Introduction

Complex analysis of large volumes of data is not unusual in MS analysis. However, signal processing approaches of the raw data in order to extract by the best of the acquisition data, has been only partly explored. In the recent years, in particular thanks to advances in mathematics, new approaches related to Compressed Sensing or to Machine Learning, allow to the development of improved sensibility and results, or even to explore new acquisition procedures.

Methods

As an example of such optimisation, the natural sparsity of MS data can be exploited through the minimization of the ℓ_1 or ℓ_2 norms. Proximal algorithms can be then developed which rapidly converge without the need to compute a (costly) gradient. Discreet pattern dictionaries and pattern matching algorithms will be presented for the analysis of isotopic patterns, and fast rank reduction operation as a method to reduce the noise, or even complement missing values.

Results

The basis of the approach will be presented on resolution enhancement of Orbitrap experiments, obtained by a fast optimization of a sparsity score function, using proximal operators. These approaches being fast enough to allow a near-real time analysis.

The same approach is then applied to FT-ICR data-sets, in particular to the 2D FT-ICR method[1]. Here, sparsity is combined to random projections to remove unwanted noise[2], enhance resolution, and missing point reconstruction[3].

Finally, we present recent work where isotopic patterns in complex are recognized in overlapped spectra, without the need of separated peaks. The use of the averaging model[4] allows to build a fast approximation of the patterns, and a search in the Fourier space allows to considerably speed-up the analysis and reduce the memory burden[5]. Promising results are already obtained with this method on peptide mixtures with several charge states intermixed, or with large proteins for which precise determination of the monoisotopic value is obtained.

Conclusions

The methods presented here rely on the analysis of the raw data, as obtained before any further processing. They have to be inserted before the usual data analysis pipeline. This displaces the chemist habits but is certainly worth the burden. This research is from a joint action with the CASC4DE spin-off company. This company works on data analysis in the biosciences, and is chosen to develop the data analysis pipeline of the EU_FT-ICR_MS European infrastructure.

Novel Aspect

New mathematics in data processing can bring you unexpected results.

References

1. Agthoven, M. A., Delsuc, M.-A., Bodenhausen, G., Rolando, C. *Anal. Bioanal. Chem.*, 405, 51–61 (2013).
2. Chiron, L., van Agthoven, M. A., Kieffer, B., Rolando, C., Delsuc, M.-A. *Proc. Nat. Acad. Sci*, 111, 1385–1390 (2014)
3. Bray, F., Bouclon, J., Chiron, L., Witt, M., Delsuc, M.-A., Rolando, C. *Anal. Chem.* 89 (17) (2017).
4. M.-W Senko, J.-P Speir, F.-W McLafferty, *Anal. Chem.* 66, 2801–2808, (1994).
5. Cherni, A., Chouzenoux, E., Delsuc, M.-A. *ICASSP conference*(2018).

FOR.02-2 ADVANCED MASS SPECTROMETRY METHODS FOR UNRAVELLING THE MECHANISMS OF ORGANOPHOSPHATE TOXICITY

Huixin Wang⁽¹⁾ - **Michael Leeming**⁽²⁾ - **Blake Cochran**⁽³⁾ - **William Donald**⁽¹⁾

⁽¹⁾ School of Chemistry, UNSW, School of Chemistry, University of New South Wales, Sydney - ⁽²⁾ School of Chemistry, The University of Melbourne, Melbourne - ⁽³⁾ School of Medical Science, University of New South Wales, Sydney

Keywords: Organophosphate; protein-target identification; bottom-up proteomics; stable isotope labelling; metabolomics;

Introduction

Organophosphates (OPs) are widely used as insecticides. Some of the most toxic organophosphates have been used in warfare causing thousands of fatalities. Acute and chronic exposure to organophosphate insecticides can cause adverse health effects. The mechanism of delayed toxicity is thought to involve off-target inhibition of serine proteases [1], although the precise molecular details remain unclear due to the lack of an analytical method for global detection of proteins that have been covalently modified by organophosphates.

Methods

A novel mass spectrometry method was developed to identify OP-adducted proteins in a non-targeted fashion. Human plasma was incubated with a mixture of isotope labelled and unlabeled OP pesticides. Proteins were extracted, digested and analyzed by LC-MS/MS to detect “twin

ions" of peptides adducted by OPs. The LC-MS/MS data was processed by a blended data analytics software, which can identify the amino acid residue of proteins modified by the labelled OPs [2-4].

Results

We found that proteins can be modified covalently by OPs by both trans-methylation and phosphorylation. Transmethylation at Glu, Asp, Cys, His and Lys residues are identified as major reactive pathways which have not been reported in the literature to date. Seven protein targets of OPs from human plasma were identified. These amino acid residue target sites were also confirmed by NMR by reacting four commonly used OP insecticides with single amino acids.

Conclusions

We discovered a new non-selective reaction pathway between OPs and human proteins. This transmethylation reaction may be responsible for the mechanism of delayed toxicity owing to organophosphate exposure. It is anticipated that this information will be useful for revealing the mechanism(s) of organophosphate toxicity.

Novel Aspect

We report the first non-targeted identification of the proteins that react relatively non-specifically with OPs. This method is more reliable and efficient compared to the traditional pull-down method [5]. Transmethylation is identified as a new reaction pathway in addition the well-known phosphorylation mechanism for OP-protein interactions.

References

1. Casida J., Quistad G., *Chem. Res. Toxicol.*, 17, 983 (2004).
2. Leeming M., Isaac, A., Pope, B., Cranswick N., Wright C., Ziogas J., O'Hair R., Donald W., *Anal. Chem.*, 87, 4104 (2015).
3. Leeming M., Donald W., O'Hair R., *Anal. Chem.* 89, 5748 (2017).
4. Hoi K., Daborn P., Battlay P., Robin C., Batterham P., O'Hair R., Donald W., *Anal. Chem.* 86, 3525 (2014).
5. Kidd D., Liu Y., Cravatt B., *Biochemistry.* 40, 4005 (2001).

FOr.02-3 ALGORITHM FOR PROCESSING COMPLEX MIXTURE DATA FROM HYPHENATED ULTRAHIGH RESOLUTION MASS SPECTROMETRY

Rémy Gavard⁽¹⁾ - **Hugh Jones**⁽²⁾ - **David Rossell**⁽³⁾ - **Simon Spencer**⁽⁴⁾ - **Mark P. Barrow**⁽²⁾

⁽¹⁾ University of Warwick, MAS CDT, Coventry - ⁽²⁾ University of Warwick, Department of Chemistry, Coventry - ⁽³⁾ Universitat Pompeu Fabra, Department of Economics and Business, Barcelona - ⁽⁴⁾ University of Warwick, Department of Statistics, Coventry

Keywords: FTICR, GC, complex mixtures, data analysis, algorithm

Introduction

The coupling of gas chromatography (GC) to a Fourier transform ion cyclotron mass spectrometer (FTICRMS) affords additional information about complex mixtures. Co-eluted components can be resolved due to the ultrahigh resolving power, also allowing extracted ion chromatograms (EICs) to be used for the observation of isomers. There is a current need for data analysis tools for processing these datasets, aiding data analysis, and extraction of key information.

Methods

A petroleum-related sample was analyzed using GC coupled with FTICR MS. An algorithm was used for pre-processing of the data. Peaks were filtered with respect to m/z and retention time

(RT) and an averaged list was generated prior to compositional assignments using the petroleomics software, Composer. The algorithm was subsequently used to incorporate the assignments into the processed data, generating a series of interactive plots.

Results

The typical workflow for this type of data is based upon manually merging data over a range of retention times, extraction of a peak list, performing the compositional assignments, visualizing the results, and repeating for each retention time range. An environmental sample from the oil sands region of Alberta, Canada, was used as a testbed for a new algorithm, which was developed using the open source R language, coupled with the array of tools offered by the Tidyverse and Shiny packages. The algorithm was used to process the dataset, monitor peaks with respect to m/z and retention time, filter peaks of inconsistent m/z , and generate an averaged peak list of all the EICs detected. The data was then visualized according to heteroatom class, carbon number, double bond equivalents (DBE), and retention time. The algorithm also affords the abilities to screen for isomeric contributions and to follow homologous series and compound classes, instead of individual components, as a function of time.

Conclusions

Using a new algorithm, a time-consuming workflow has been replaced with one taking only a few minutes. It has become possible to follow both individual components and compound classes as a function of retention time, without the necessity for averaging multiple time windows and the accompanying loss of information. The algorithm is not only suitable for application to GC data, but also data produced using other separation techniques.

Novel Aspect

New algorithm developed in R to process and visualize complex mixture data, acquired using hyphenated ultra high resolution mass spectrometry.

FOr.02-4 **LC-MS1 METHOD FOR PROTEIN IDENTIFICATION AND QUANTITATION WITHOUT USING TANDEM MASS SPECTRA**

Mark Ivanov - Irina Tarasova - Elizaveta Solovyeva - Mikhail Gorshkov

V.L. Talrose Institute for Energy Problems of Chemical Physics Russian Academy of Sciences, Moscow

Keywords: proteomics, protein identification, LFQ quantitation

Introduction

A major bottleneck of modern shotgun proteomics is limitation of fragment mass spectra acquisition speed. We present a novel workflow for protein identification based on LC-MS1 analysis with no using of tandem mass-spectra [1]. The proposed method has no drawbacks of accurate mass tags (AMT) or classic peptide-mass fingerprinting approaches.

Methods

Proposed approach combines prediction of peptide retention times and accurate peptide mass measurements obtained using high resolution mass spectrometers. Dinosaur software [2] is used for peak picking, ELUDE algorithm [3] is used for retention time prediction and Diffacto algorithm [4] is used for quantitation analysis.

Results

Here we present identification and quantitation efficiency of developed workflow. The MS1-based approach outperforms traditional MS/MS search in both number of identified proteins groups and sequence coverage for the HeLa cell line analyzed on Orbitrap Velos with 25 min LC

gradient. We have shown that LC-MS1 based approach does not depend on the length of chromatographic gradient and can potentially replace traditional MS/MS approach in high throughput proteomics applications. Also, we have shown that MS1-based workflow successfully found changes in concentration of all 6 proteins from iPRG2015 data set without any false positives. For the blood plasma samples, the MS1-based approach reported concentration changes in 4 proteins which are known FDA-approved disease markers.

Conclusions

The ability of using the developed method for short gradients with better efficiency compared with MS/MS can potentially be an interesting feature for the applications requiring high throughput analyses of clinical samples using relatively simple and inexpensive LS-MS setups.

Novel Aspect

A novel scheme for successful protein identification in complex mixtures using short LC gradients and no peptide fragmentation.

References

1. Ivanov M.V. et al., *J. Proteome Res.*, 16, 3989-3999 (2017)
2. Teleanu J. et al., *J. Proteome Res.*, 15, 2143-2151 (2016)
3. Moruz L. et al., *J. Proteome Res.*, 9, 5209-5216 (2010)
4. Zhang B. et al., *Mol Cell Proteomics*, 16, 936-948 (2017)

FOr.02-5 PHASE-CONSTRAINED SPECTRUM DECONVOLUTION–REAL TIME SUPER-RESOLUTION PROCESSING OF FTMS DATA

Konstantin Aizikov - Dmitry Grinfeld - Arne Kreuzmann - Daniel Mourad - Eugen Damoc - Tabiwang N. Arrey - Oliver Lange - Alexander Makarov

ThermoFisher Scientific, LSMS, Bremen

Keywords: FTMS, Orbitrap, super-resolution, high-throughput, TMT

Introduction

The **Phase-Constrained Spectrum Deconvolution Method** (Φ SDM) [1] is a novel signal analysis approach which enables qualitatively faster FTMS scan rates while preserving the levels of mass resolution and accuracy. This presentation discusses real-time Φ SDM processing of complex experimental spectra that provide the increased depth of analyte coverage.

Methods

The MS experiments were conducted on the research grade Q Exactive™ family of Orbitrap™ mass spectrometers. The LC separation was done on EASY-Spray™ PepMap™ RSLC C18 column using UltiMate™ 3000 NCS system. Intact proteins were infused directly. An auxiliary computer equipped with four Titan X Nvidia™ GPU cards was used for the full mass range Φ SDM calculations. Spectral windows comprising the TMT labels were processed directly on the instrument's computer.

Results

Back to back LCMS runs of 11-plex TMT labeled HeLa digests demonstrated that the Φ SDM facilitates a significant increase in the spectral acquisition rate on a research grade Q Exactive HF which results in 1.74x more MS/MS scans and allowed the 1.55x PSM increase. The numbers of quantified peptides and proteins is boosted by factors of 1.47 and 1.32, respectively. The Q Exactive HF-X data shows the 2.39x MS/MS scan rates and 1.99x PSM's. The peptide and protein quantification improvement was demonstrated as 1.82x and 1.62x, correspondingly.

Real-time Φ SDM processing of intact protein spectra isotopically resolves the entire charge state envelope of Ubiquitin (8.6 kDa) with MS scan rate in excess of 30 Hz, with the reported deconvolved nominal mass accuracy of \sim 1ppm. For the carbonic anhydrase (29 kDa), the baseline isotopic resolution of the charge state envelope is achievable at scan rates of \sim 7.5 Hz, with the reported nominal mass accuracy of \sim 2ppm.

Conclusions

Φ SDM has been used for real-time LC experiments. Its ability to provide the super-FT resolution allowed higher MS/MS rates and, consequently, boosted up the identification and quantification. Protein isotopic envelopes were completely resolved with Φ SDM even at high repetition rates, which allowed reliable charge state identification and accurate reporting of the nominal mass.

Novel Aspect

Φ SDM, a new FTMS algorithm, was for the first time implemented for real-time data processing. Quantitation improvement and superior mass resolution are demonstrated.

References

1. G. Grinfeld *et al.*, *Anal. Chem.*, **2017**, *89*, 1202–1211
2. A. Thompson *et al.*, *Anal. Chem.*, **2003**, *75*, 1895–1904.

FOr.03 - Cultural heritage and archaeology

Chairs: Klaas Jan van den Berg, Giorgio Mellerio

FOr.03-1 ***Keynote: ART AND CULTURAL HERITAGE NATURAL POLYMERS BY BOTTOM UP AND TOP DOWN APPROACHES***

Caroline Tokarski

University of Lille

Keywords: Cultural Heritage, Top down proteomics, Lipidomics, Metabolomics, Crosslinkings

Introduction

Natural organic polymers represent the main organic components of the most of the historic samples. However Cultural Heritage material is challenging to study[1]: (i) The sample size typically obtained from a museum- or archaeological artefact is on the order of tens of micrograms. (ii) Most pieces are made of more than one material, each containing many organic/inorganic components. (iii) The material is constantly evolving due to ageing or external conditions.

Methods

Top down experiments are used for protein studies; particular focus is given to sample preparation fitted to handle trace compounds. Soft depolymerization and HRMS analysis of multimers are used for lipids; Kendrick and Van Krevelen plots help to decipher the high complexity of spectra. Polysaccharide analysis combines filter aided sample preparation and nanoLC-HRMS analysis of fingerprints; methods are fitted both to linear and highly branched structures

Results

This conference will show how organic chemistry combined with HRMS (orbitrap, FT-ICR) address the current challenges for analysis of trace amounts of Cultural Heritage samples. For

example structural elucidation of biopolymers of unknown structures and/or their chemical modifications will be shown (e.g. historic art paintings and watercolors [2]). Focus on protein crosslinkings inside paintings will be done. Another example is the combination of soft depolymerization experiments and HRMS to unravel the 3D networks formed by insoluble lipidic films (e.g. oil-paintings). The very first application of top down approach to the study of protein remains in archaeological ceramics will show impressive results both in terms of protein sequence coverage (up to 100% coverage) and identified sequences that show particular degradation patterns. The conference will be illustrated by the study of various outstanding samples from art but also from archaeological and paleontological sciences from the most famous museums and institutes of the world.

Conclusions

Applied to the study of artworks, the chemical decoding of the organic molecular network has a major impact on the understanding of painting formulation, biopolymers degradation and it helps in the survey of preservation conditions. Considering archaeological samples, the information obtained by the proposed methods provide new historical insights; for example it points out ancient methods for food processing.

Novel Aspect

Deciphering Art and Cultural Heritage natural polymers using unexplored developments in chemistry and top down approaches to chemically decode the organic molecular network

References

1. Dallongeville S., Garnier N., Rolando C., Tokarski C., *Chemical Reviews*, 116, 2-79 (2016).
2. Granzotto C., Arslanoglu J., Rolando C., Tokarski C. *Scientific Reports* 7, 44538 (2017).

FOR.03-2 FAST ANALYSIS OF MODERN OIL PAINTINGS AND DYES IN TEXTILES BY SURFACE ACOUSTIC WAVE NEBULIZATION-MASS SPECTROMETRY (SAWN-MS)

Alina Astefanei⁽¹⁾ - ***Maarten van Bommel***⁽²⁾ - ***Klaas Jan van den Berg***⁽³⁾ - ***Garry Corthals***⁽¹⁾

⁽¹⁾ University of Amsterdam, HIMS, Analytical Chemistry, Amsterdam - ⁽²⁾ Cultural Heritage Agency of the Netherlands, and HIMS, Analytical Chemistry, University of Amsterdam, Amsterdam - ⁽³⁾ Cultural Heritage Agency of the Netherlands (RCE), and University of Amsterdam, Amsterdam

Keywords: SAWN-MS, oil paintings, dyes, textiles

Introduction

GC-MS and LC-PDA are the methods of choice to study the binding media of paintings and dyes in textiles [1,2]. Relatively large samples are needed and long sample preparation and analysis time. In GC-MS, limited information is obtained due to the hydrolysis and derivatization steps, the high temperatures and ionization energies used.

SAWN-MS maintains the integrity of fragile molecules due to mild ionization, greatly reduces the sample sizes and the analysis time.

Methods

The SAWN device was placed in front of the API inlet of a triple-TOF 5600+ MS (Sciex) and operated through a connection with a frequency and power controller. 1 μ L of liquid sample was loaded on the chip and the SAW regulated by applying a power to the electrodes (\approx 5 W). This step was repeated if needed (up to 1 min) and the data accumulated into a single file. The interface heater temperature was 150 $^{\circ}$ C, the inlet/outlet and curtain gasses were 0 psi and 10 psi respectively.

Results

Herein, we describe the SAWN-MS technique and optimised experimental procedures for the analysis of modern oil paints and dyes in textiles for conservation studies.

SAWN-MS analysis revealed a rich representation of molecular ions that were used for the identification of mixtures of colourants present in dyed textiles [3] providing information on the colouring procedure and the state of degradation.

The technique was also applied to study levels of oxidation and hydrolysis of 20th century oil paints, by monitoring the extracted free fatty acids and diacids, both as glycerides and in the free form. The formation of diacid functionalities play a significant role in the development of water sensitivity and paint instability/failure [4]. A correlation with the pigment type and additives used in paint formulations was obtained. SAWN-MS also enabled the detection of short chain fatty acids, which cannot be easily detected by GC-MS.

Conclusions

SAWN-MS is an excellent analytical tool for conservation studies on oil paintings and textiles. It offers unique advantages in terms of ease of implementation, simplifies the sample extraction procedures, significantly reduces the required sample sizes and overall analysis times are typically few minutes. Moreover, due to the new mechanisms and efficiency of ionization it provides additional molecular information compared to other ionization techniques.

Novel Aspect

Novel ionization technique for fast MS analysis of organic components in paintings and dyes in textiles.

References

1. Bonaduce, I., Ribechini, E., Modugno, F., Colombini, M. P., In *Analytical Chemistry for Cultural Heritage*; Mazzeo, R., Ed. Springer International Publishing: Cham, 291–327 (2017).
2. Valianou, L., Karapanagiotis, I., Chryssoulakis, Y., *Anal. Bioanal. Chem.*, 395 (7), 2175–2189 (2009).
3. Astefanei, A., van Bommel, M., Corthals, G. L., *J. Am. Soc. Mass Spectrom.*, 28 (10), 2108–2116 (2017).
4. Banti, D., La Nasa, J., Lluveras Tenorio, A, Modugno, F., van den Berg, K. J., Lee, J., Ormsby, B., Burnstock, A., Bonaduce, I. *RSC Adv.*, 8, 6001, (2018).

FOr.03-3 **MOLECULAR SIGNATURES OF PROTEIN AGING IN ANCIENT HUMAN BONES FROM THE ERUPTION AREA OF VESUVIUS**

Georgia Ntasi⁽¹⁾ - **Piero Pucci**⁽¹⁾ - **Gennaro Marino**⁽²⁾ - **Leila Birolo**⁽¹⁾

⁽¹⁾ *University of Naples Federico II, Chemical Sciences, Naples* - ⁽²⁾ *University of Naples, Suor Orsola Benincasa, Naples*

Keywords: Paleoproteomics, Mass Spectrometry, Cultural heritage, Protein Degradation, Ancient proteins

Introduction

In search of markers for protein aging, we use the singular case of collagen from human bones from the excavations sites in Pompei, Herculaneum, Oplonti, S. Paolo Belsito, and Scalandrone Bay as control reference site. Deamidation of glutamine and asparagine is considered as aging molecular marker [1]. Moreover, other chemical modifications that are a likely signature of diagenetically derived alterations were considered.

Methods

Specific proteomic protocols and different sample preparation were adapted to improve protein extraction and sequence coverage yields. The analyses were performed with the use of LC-MS/MS. The exploitation of bioinformatic tools such as the error tolerant search in the Mascot suite [2] or the open mass modification utility in the Protein Prospector database search [3] were used for data processing in search of unexpected chemical modifications.

Results

Ancient human bones (100) have been analyzed by proteomic approach following an adapted protocol based on the already developed one reported in [4]. From the raw LC-MS/MS data and the use of Mascot platform, COL1 α 1 and COL1 α 2 were detected with a protein sequence coverage up to 80% in all the samples. All the deamidation sites were listed and their occurrence and deamidation state were statistically evaluated. The statistical analysis is suggestive of the most detected positions, that were considered the best candidates for the quantification of the deamidation level as protein's aging biomarkers in ancient human bones. Similar approach was used to evaluate also other chemical modifications.

Conclusions

The analysis of the deamidation in the different samples suggests that the deamidation level in the samples from Pompei is low in respect to an extensive deamidation observed in the samples from Herculaneum and Scalandrone (our control reference). Our preliminary results seem to suggest that the low deamidation level could be linked to the peculiar burial conditions and to the "lower" temperature experienced by Pompei inhabitants in comparison to Herculaneum ones.

Novel Aspect

This set of samples constitutes an extraordinary test case to evaluate deamidation and other chemical modifications as molecular signatures of protein aging.

References

1. Leo G, Cartechini L, Pucci P, Sgamellotti A, Marino G, Birolo L. Proteomic strategies for the identification of proteinaceous binders in paintings. *Anal. Bioanal. Chem.* 2009;395(7):2269-2280. doi:10.1007/s00216-009-3185-y.
2. Cappellini E, Jensen LJ, Szklarczyk D, et al. Proteomic Analysis of a Pleistocene Mammoth Femur Reveals More than One Hundred Ancient Bone Proteins. *J. Proteome Res.* 2012;11(2):917-926. doi:10.1021/pr200721u.
3. Hill RC, Wither MJ, Nemkov T, et al. Preserved Proteins from Extinct Bison latifrons Identified by Tandem Mass Spectrometry; Hydroxylysine Glycosides are a Common Feature of Ancient Collagen. *Mol. Cell. Proteomics.* 2015;14(7):1946-1958. doi:10.1074/mcp.M114.047787.
4. Vinciguerra R, De Chiaro A, Pucci P, Marino G, Birolo L. Proteomic strategies for cultural heritage: From bones to paintings. *Microchem. J.* 2016;126:341-348. doi:10.1016/j.microc.2015.12.024.

For.03-4 CHARACTERIZATION OF ISOBARIC SPECIES AND PHOTO-DEGRADATION PRODUCTS OF TRIARYLMETHANE DYES WITHIN TEXTILES

Francesca Sabatini⁽¹⁾ - **Chiara Bracini**⁽¹⁾ - **Anita Quye**⁽²⁾ - **Ilaria Degano**⁽¹⁾ - **Maria Perla Colombini**⁽¹⁾

⁽¹⁾ University of Pisa, Department of Chemistry and Industrial Chemistry, Via Giuseppe Moruzzi 3, 56124 Pisa, Italy - ⁽²⁾ Centre for Textile Conservation and Technical Art History, School of Cultural and Creative Arts, University of Glasgow, Glasgow G12 8QH, United Kingdom
*f.sabatini4@gmail.com

Keywords: LC-ESI-Q-ToF, Triarylmethane dyes, accelerated ageing, photo-oxidation, degradation pathways

Introduction

Since mid-19th century, triarylmethane dyes were widely used for ink production and dyeing textiles [1]. These dyes, characterized by bright hues, are one of the first classes of synthetic dyes produced and marketed [2]. Nevertheless, their analysis is highly complicated due to their poor light fastness and their complex composition constituted by homologous species, often differing only for the presence or position of the same substituent on the aromatic rings.

Methods

Liquid Chromatography coupled with tandem high-resolution mass spectrometry (LC-ESI-Q-ToF) or Diode Array and Fluorescence detectors (HPLC-DAD-FD) were used to separate and identify the several miscellaneous components contained in standard triarylmethane dyes, unaged and aged wool yarns and historical samples of dyed specimens. The samples were pre-treated using an oxalic acid solution to obtain the highest extraction yields for all the components.

Results

The analysis of standards and reference materials by LC-ESI-Q-ToF allowed us to identify specific fragmentation pathways for several triarylmethane dyes such as Fuchsine, Crystal/Methyl Violets, Methyl Blue and Diamond Green. Both positive and negative modes and collision energies were tested to optimize the MS conditions for each dye and to obtain complementary information on fragmentation pathways.

The analysis of wool yarns dyed with triarylmethane compounds, artificially aged for different time intervals, allowed us to investigate and highlight different degradation pathways, through the identification of kinetic trends of different compounds in relation to the accelerated aging times. HPLC-DAD-FD was fundamental to provide further information on the minor compounds. Analysis of selected historical samples from 19th C dyeing manuals from Glasgow's 'Dye-versity' project [3] highlighted the compositional complexity of early synthetic dyes, and revealed several isomers and different profiles depending on desired hues.

Conclusions

The proposed analytical method, integrating HPLC-ESI-Q-ToF and HPLC-DAD-FD analyses, allowed the full characterization of the references and historical samples under study and to hypothesize the degradation pathway of several triarylmethane compounds. The detection of specific fragmentation patterns enabled distinction between different classes of compounds, while that of specific fragment ions distinguished isobaric species belonging to the same series.

Novel Aspect

We used for the first time a solid silica core analytical column coupled with HRMS to fully characterize minor and degradation components in reference and historical samples.

References

1. Confortin D., Neevel H., Brustolon M., Franco L., Kettelarij A.J., Williams M.R., van Bommel M.R., *Journal of Physics: Conference Series*, 231, 1-9 (2010)
2. Fyfe I., *Arts & Reviews: Modern World, Modern Art*, Bluesci-Cambridge University Magazine, Vol. 19, p. 26-27, 2010.
3. Quye, A., Han. J. 'Disguises, surprises and identity crises for 19th century early aniline dyes'. Paper presented at the 35th Annual meeting of Dyes in History and Archaeology, 5-8 October 2015, Pisa. <https://www.gla.ac.uk/schools/cca/research/arthistoryresearch/projectsandnetworks/dye-versity/>

For.03-5 WINE PRODUCTION AND CONSUMPTION IN PRE ROMAN TIMES: RESIDUE ANALYSIS OF A PRODUCTION STRUCTURE AN CERAMIC VESSELS FROM THE WESTERN MEDITERRANEAN

Alessandra Pecci⁽¹⁾ - ***Andrea Zifferero***⁽²⁾ - ***Gianluca Giorgi***⁽³⁾ - ***Simona Mileto***⁽¹⁾ - ***Rafael Jornet***⁽¹⁾ - ***Carme Belarte Franco***⁽⁴⁾ - ***Anna De Palmas***⁽⁵⁾ - ***Massimo Vidale***⁽⁶⁾

⁽¹⁾ ERAAUB, Departament Història i Arqueologia, Universitat de Barcelona, Spain - ⁽²⁾ Department of History and Cultural Heritage, University of Siena, Italy - ⁽³⁾ Department of Biotechnology, Chemistry and Pharmacy, University of Siena, Italy - ⁽⁴⁾ CSIC, Spain - ⁽⁵⁾ DUMAS, Università di Sassari - ⁽⁶⁾ Department of Cultural Heritage, University of Padua (Italy)

Keywords: organic residue analysis, tartaric acid, wine, gas chromatography-mass spectrometry

Introduction

Wine was an important product in the Western Mediterranean. Organic residue analysis was carried out on samples from pottery and a production structure in four pre - Roman archaeological sites in Italy and Spain, aimed at identifying the presence of wine. At most of the sites the production and consumption of wine is compatible with other archaeological evidence (mainly botanical evidence).

Methods

Samples were recovered both from the archaeological materials and the earth of the stratigraphic unit to which it belonged, in order to check possible contamination.

The samples were mechanically cleaned in the laboratory and powdered.

Samples were extracted following the extraction method proposed by Pecci et al. [1], derivatized and analysed by gas chromatography coupled with mass spectrometry.

Results

The analyses indicate the presence of tartaric acid in some of the analysed pottery samples from the different sites in Italy and in a sample from an Iberian production structure at Coll del Moro (Catalunya, North East Iberian Peninsula).

Although other fruits, for example tamarind, may contain tartaric acid [2], on the basis of the residues identified [1,3,4], and the presence of grape pips or pollen in the area it is possible to suggest that wine or its derivatives was in contact with the sampled materials.

The fact that no tartaric acid was identified in control samples, nor in some of the archaeological samples indicates that it does not derive from post depositional contamination.

In some of the samples also other acids such as succinic, fumaric, isocitric, malic and maleic acids, are present. Although they are not exclusive markers of wine, they are present in the beverage.

Conclusions

The data in our possess allow to confirm the hypothesis on the production of wine in the Iberian structure at Coll del Moro where abundant grape pips had been identified and the storage and/or consumption of wine or its derivatives at Etruscan and Bronze Age sites in Italy.

Novel Aspect

Investigation on wine production and consumption, not only in ceramic vessels, but also in a pre-Roman production structure.

References

1. Pecci A., Giorgi G., Salvini L., Cau M. A., Journal of Archaeological Science, 40, 109-115(2013).
2. Barnard H., Dooley A.N., Areshian G., Gasparyan B., Faull K.F., Journal of Archaeological Science, 38, 977-984 (2010).

3. Mc Govern P., Jalabadze M., Batiuk S., Callahan M.P., Smith K., Hall G., Kvavadze E., Maghradze D., Rusishvili N., Bouby L., Failla O., Cola L., Mariani L., Boaretto E., Bacilieri R., This P., Wales N., Lordkipanidze L., PNAS, doi/10.1073/pnas.1714728114 (2017).
4. Garnier N., Valamoti S., Journal of Archaeological Science, 76, 195-206 (2016).

FOr.04 - Atomic mass spectrometry: metallomics and isotope ratiomass spectrometry

Chairs: Ryszard Lobinski, Federica Camin

FOr.04-1 **Keynote: NEW ION SOURCES FOR METALLOMIC ANALYSIS AND ISOTOPE-RATIO MASS SPECTROMETRY**

Gary M. Hieftje

Department of Chemistry - A150 Chemistry Building, Indiana University, Bloomington

Keywords: metallomics, speciation, plasma source MS, Isotope ratio, instrumentation

Introduction

Most ion sources for mass spectrometry are tailored for specific species such as molecules, atoms, or biomolecules, requiring instruments to be designed and constructed with limited flexibility. What would be desirable is a source that can ionize a broad range of species and sample types, including solids, liquids, or vapors. Such a source, termed the Solution-Cathode Glow Discharge (SCGD) is now being evaluated and will be described and evaluated

Methods

The Solution-Cathode Glow Discharge (SCGD) is a relatively low-power (~70W) ionization source that operates in the open atmosphere. Sample solution flows from a small glass capillary and serves as the cathode of a discharge that strikes to a metallic anode pin just above it. The electric field at the cathode surface is high enough that electrical spraying takes place. Thermal fragmentation or even complete atomization of incoming analyte species is also possible.

Results

A broad range of sample types has been investigated with the SCGD. Metallic species in solution can be either fully atomized or partially fragmented, yielding alternatively atomic mass spectra or spectra reflecting original solution composition. Thus, the SCGD is useful for either elemental or chemical-speciation analysis. Further, sensitivity is sufficient that isotope ratios can be measured at precision levels dictated by the IAEA International Target Values for natural uranium (0.28% rsd), with an integration time of only 150 s. For somewhat labile molecules, fragmentation can be controlled merely by adjusting the operating current and voltage of the SCGD. For triethyl phosphate, for example, the survival yield can be varied from at least 25-80%. Under conditions that favor electrical spraying, the SCGD can be fed directly with a peptide-containing solution and yields fragmentation mass spectra that are unlike those from other ionization sources and are surprisingly easy to interpret.

Conclusions

Although originally developed as an emission source for elemental analysis, the SCGD has proven to be an extremely versatile ionization source for mass spectrometry. It is simple to construct, operates on a low-power supply, requires no gas supply, and is simple to set up and employ. Most work to date has been with an Orbitrap® spectrometer; ongoing work will employ alternative mass spectrometers, some suitable for field use.

Novel Aspect

A new, highly versatile but simple ion source has been developed and is being tested. It is useful for samples ranging from biomolecules to atoms.

FOr.04-2 ELEMENT ANALYSIS WITH THERMALIZED LASER ABLATION IONS FROM SOLID SAMPLES

Lorenzo Querci - Bodo Hattendorf - Detlef Günther

ETH Zürich, Laboratory of Inorganic Chemistry, Zurich

Keywords:

Laser ablation/ionization mass spectrometry, collisional cooling, ion funnel, fs-laser, CD-nozzle

Introduction

Laser ablation/ionization mass spectrometry is becoming again an attractive analytical technique for element composition determination of solids [1- 4]. It offers several advantages: little sample preparation, high sample throughput, high spatial resolution and the possibility to analyze any kind of solid sample. We present a novel strategy using fs-LA at elevated pressure in combination with CD nozzle and Rf-only ion funnel for element imaging in solids.

Methods

The laser ablation ion funnel ion source, described in this work, employs an Rf only ion funnel in combination with a convergent-divergent nozzle for collisional cooling of laser generated ions and their separation from the buffer gas. This configuration achieves a sufficient pressure drop between ion source and mass spectrometer while allowing for multiple collisions and consequent focusing. The "LAFU" source was coupled to a TOF-MS and different laser were employed.

Results

The first results obtained with the coupling of ns-LAFU (3 ns, 532 nm) source with time of flight mass spectrometry (LAFUMA), showed that also low m/Q ions (^{12}C , ^{27}Al) can be successfully focused in the Rf-only ion funnel. The ion signal intensities could be increased by up to 300 times in comparison with ablation in high vacuum. Still the ion transmission through the ion funnel was found to have a pronounced dependence on m/Q , as a result of the pseudo-stationary fields created near the ion funnel electrodes. Compared to the previously used ns laser ablation, the fs laser (150 fs, 400 nm) ablation exhibits significantly smaller mass ablation rates and generally more comparable sensitivities. Transient evolution of the pulsed ion signals was longer than for ns LA regime and less affected by buffer gas pressure., indicating that a fraction of high energetic ions could not be thermalized completely in this case.

Conclusions

An ion funnel has been successfully combined with laser ablation for ion generation at elevated pressure, achieving a significant enhancement of ion transmission. The combination of the ion funnel with a CD nozzle for differential pumping shows promising features for the further development of laser ionization/ablation mass spectrometry in element and isotope analyses of solids.

Novel Aspect

Collisional cooling of laser generated ions at elevated pressure in combination with an Rf only ion funnel for fast element image acquisition was successfully realized.

References

1. R. E. Honig, *Appl. Phys. Lett.*, 3, 8 (1963).
2. N. C. Fenner and N. R. Daly, *Rev. Sci. Instrum.*, 37, 1068–1070 (1966).

3. A. Riedo, M. Neuland, S. Meyer, M. Tulej and P. Wurz, *J. Anal. At. Spectrom.*, 28, 1256. (2013).
4. M. He, B. Li, B. Zhang, W. Hang and B. Huang, *J. Anal. At. Spectrom.*, 27, 2–7 (2012).

FOr.04-3 HIGH-RESOLUTION AND TANDEM MASS SPECTROMETRY TECHNIQUES CONTRIBUTE TO THE DISCOVERY OF NEW ARSENIC METABOLITES OF TOXICOLOGICAL IMPORTANCE

Qingqing Liu⁽¹⁾ - Xiufen Lu⁽¹⁾ - Hanyong Peng⁽¹⁾ - Elaine Leslie⁽²⁾ - Martin Zuidhof⁽³⁾ - Xing-Fang Li⁽¹⁾ - X. Chris Le⁽¹⁾

⁽¹⁾ University of Alberta, Analytical and Environmental Toxicology, Edmonton - ⁽²⁾ University of Alberta, Physiology, Edmonton - ⁽³⁾ University of Alberta, Agricultural, Life and Environmental Sciences, Edmonton

Keywords: Arsenic speciation; protein binding; toxicology; biotransformation; LC/MS/ICPMS

Introduction

Nearly 100 arsenic species, with diverse toxicities (1,2), are present in the environment and in biological systems. There is a tremendous analytical challenge in the identification and quantification of trace amounts of individual arsenic species present in complex biological sample matrix. This presentation will highlight our recent development of mass spectrometry techniques for studying arsenic species, including the oxygenated and thiolated, trivalent and pentavalent, and inorganic and organic arsenic species.

Methods

Hybrid quadrupole time-of-flight, linear ion trap, and triple quadrupole mass spectrometers were used to conduct high-resolution and tandem mass spectrometry experiments. Effluent from high performance liquid chromatography (HPLC) separation was split 80% to inductively coupled plasma mass spectrometry (ICPMS) and 20% to electrospray ionization tandem mass spectrometry (ESI-MS/MS). Human primary hepatocytes and hepatocellular carcinoma HepG2 cells were treated with 0, 20, or 100 μM Roxarsone (3-nitro-4-hydroxyphenylarsonic acid).

Results

The combined techniques enabled identification and quantitation of individual arsenic species present in a variety of environmental and biological samples. Roxarsone was metabolized to more than 10 arsenic species in human hepatic cells. Several new metabolites were identified (3,4). The 24-hour IC_{50} values of thiolated Roxarsone for A549 lung cancer cells and T24 bladder cancer cells were 380 ± 80 and 42 ± 10 μM , respectively. This thiolated arsenic metabolite was more toxic than Roxarsone, whose 24-hour IC_{50} values for A549 and T24 cells were 9300 ± 1600 and 6800 ± 740 μM , respectively.

Conclusions

We have characterized the metabolism of Roxarsone in human liver cells and successfully identified the new phenylarsenical metabolites. The identification and toxicological studies of the new arsenic metabolites are useful for understanding the fate of arsenic species and assessing the potential impact of human exposure to Roxarsone

Novel Aspect

The complementary techniques enabled the discovery of new arsenic metabolites that have important toxicological implications. The study improved mechanistic understanding of arsenic biotransformation.

References

1. Moe B., Peng H.Y., Lu X.F., Chen B.W., Chen L.W.L., Gabos S., Le X.C.J. *Environ. Sci. Technol.* 49, 113-124 (2016).

2. Liu Q., Lu X., Peng H., Popowich A., Tao J., Uppal J.S., Yan X., Boe D., Le X.C. *Trends Anal. Chem.* DOI: 10.1016/j.trac.2017.10.006 (2018).
3. Peng H., Hu B., Liu Q., Li J., Li X.-F., Zhang H.Q., Le X.C. *Angew. Chem. Int. Ed.* 56, 6773-6777 (2017).
4. Liu Q., Leslie E., Moe B., Zhang H., Douglas D., Kneteman N., Le X.C. *Environ. Sci. Technol.* 52, 1386-1392 (2018).

For.04-4 ISOTOPIC TECHNIQUES FOR STUDYING THE SOURCES AND PROCESSES OF POLLUTANTS IN THE ENVIRONMENT

Nives Ogrinc

Jožef Stefan Institute Department of Environmental Sciences Ljubljana

Keywords: isotopic techniques, sources, mercury, organic pollutants, environment

Introduction

The use of isotopic techniques is in the midst of a remarkable period of innovation and discovery; the last decade has seen the emergence of stable isotopes of metals and new capabilities for measurements of organic compounds. These advances hold the potential to reveal new insights into the nature, sources and transformation processes of pollutants in the environment. The selected examples will include: polycyclic aromatic hydrocarbons (PAH), atrazine and mercury.

Methods

The sources of PAHs were determined in lake sediments using isotopic composition of carbon, while isotopic composition of nitrogen and carbon in atrazine was used to characterize degradation processes. These measurements were performed by gas chromatography coupled to isotope ratio mass spectrometry through combustion unit (GC-C-IRMS). Stable isotopes of Hg determined by MC-ICP MS were applied to determine the Hg sources in Mediterranean deep-sea sediments.

Results

The dominant PAHs signatures identified in sediments of Lake Bled were mainly attributed to a coal/wood burning source, but PAH from carsoot could also contribute. Retene (Re) and Perylene (Per) exhibit profiles and isotopic composition that are distinct from those of other PAHs suggesting their natural origin. However $\delta^{13}\text{C}$ values of Re determined in anoxic part linked Re to both natural and pyrolytic origin [1]. Stable isotope analysis of atrazine brings forward two additional lines of evidence to assess the degradation state of atrazine and desethylatrazine in natural samples: (i) to distinguish different degradation pathways of atrazine and (ii) to detect further degradation of DEA[2]. Both mass-dependent and mass-independent fractionation processes influence the isotopic composition of Hg in the Mediterranean Sea. Positive $\Delta^{199}\text{Hg}$ values are likely the result of enhanced Hg^{2+} photoreduction in the Mediterranean water column, while mass-dependent fractionation decreases $\delta^{202}\text{Hg}$ values due to kinetic isotope fractionation during deposition and mobilization. At least three primary Hg sources of atmospheric deposition could be defined: urban, industrial and global precipitation-derived.

Conclusions

The results indicate that together molecular and isotope composition of PAHs are useful tools for obtaining more detailed identification about sources of PAH inputs into lacustrine sedimentary environments, while trends in multielement isotope data of atrazine may, decipher different degradation pathways. Stable isotopes of Hg provide new insights into the sources and geochemistry of Hg in deep-sea sediments.

Novel Aspect

Isotopic techniques such as GC-C-IRMS and MC-ICP-MS open new avenues for studying sources, fate and impact of contaminants in the environment.

References

1. Gams Petrišič, M., Muri, G., Ogrinc, N., *Environ. Sci. & Technol.* 2013, 47, 1280-1286.
2. Meyer, A.H., Elsner, M., *Environ. Sci. & Technol.* 2013, 47, 6884–6891.

FOr.04-5 **POTENTIAL USE OF THE STABLE ISOTOPE RATIOS OF BIOELEMENTS AND ELEMENTAL COMPOSITION TO TRACE THE ORIGIN OF DAIRY PRODUCTS**

Luana Bontempo⁽¹⁾ - **Roberto Larcher**⁽²⁾ - **Federica Camin**⁽¹⁾

⁽¹⁾ Department of Food Quality and Nutrition, Research and Innovation Centre, Fondazione Edmund Mach (FEM), Via E. Mach 1, 38010 San Michele all'Adige (TN), Italy – ⁽²⁾ Technology Transfer Centre, Fondazione Edmund Mach, via E. Mach 1, 38010 San Michele all'Adige, Italy

Keywords: stable isotope ratios, elemental composition, traceability, dairy products

Introduction

There is an increasing demand for reliable analytical methods to verify the authenticity of the food we eat. This is also important for PDO dairy products, which can command significantly higher prices than their non-PDO competitors. Stable isotope ratios and elemental composition have been used for discrimination and can be used to verify the authenticity of marketed products. An overview of the use of this approach to trace dairy products is presented here.

Methods

Isotope Ratio Mass Spectrometry and Inductively Coupled Plasma Mass Spectrometry methods were developed to characterise and determine the authenticity of different types of dairy products. In particular, the stable isotope ratios (SIRs) of five bioelements- hydrogen, carbon, nitrogen, oxygen and sulphur - were determined in defatted dairy products, whereas elemental composition was determined in the raw product after acid microwave digestion.

Results

Geo-climatic and pedo-geological factors affect SIRs in nature and the isotopic variations are ultimately incorporated into animal tissue through eating, drinking, breathing and exchanges with the environment, and are memorised in the resulting dairy products. Similarly, the elemental composition of dairy products is affected by different factors related to the place of origin (e.g. geology and soil characteristics) and animal diet, but also depends closely on the conditions in which such products are made and processed. The combination of these two approaches has led to the development of models able to trace the origin of high premium products. The developed method was applied to test the traceability of PDO cheeses: Italian alpine cheeses (e.g. Fontina, Montasio and Toma) [1, 2], as well as buffalo Mozzarella from the Campania region, Grana Padano and Parmigiano Reggiano [3, 4, 5].

Conclusions

This work considers a robust and effective method that can be used to guarantee the authenticity of dairy products to consumers, focusing particularly on cheese. The method is based on determination of stable isotope ratios and the elemental composition of products. Specifically, in the case of PDO Grana Padano and Parmigiano Reggiano cheese, the model has been officially adopted as a reference method for assessing the authenticity of grated and shredded products on the market. Furthermore, this method was recently officially validated and adopted by UNI (Italian Standardisation Body, UNI Regulation 11692:2017).

Novel Aspect

The method applied and described has been officially validated and can be used before the courts in legal cases.

References

1. Bontempo L., Larcher R., Camin F., Hoelzl S., Rossmann A., Nicolini G., International Dairy Journal, 21, 441-446 (2011).
2. Bontempo L., Lombardi G., Paoletti R., Ziller L., Camin F., International Dairy Journal, 23, 99-104 (2012).
3. Camin F., Wehrens R., Bertoldi D., Bontempo L., Ziller L., Perini M., Nicolini G., Nocetti M., Larcher R., Analytica Chimica Acta, 711, 54-59 (2012).
4. Camin F., Bertoldi D., Santato A., Bontempo L., Perini M., Ziller L., Stroppa A., Larcher R., Rapid Communications in Mass Spectrometry, 29, 415-423 (2015).
5. EU Regulation 584/2011, Official Journal of the European Union, 160, 65-70 (2011).

FOr.05 - Nutraceuticals directions

Chairs: Roberto Romero Gonzalez, Roberta Galarini

FOr.05-1 **Keynote:** «NON SI VIVE DI SOLO PANE». DART AND AMBIENT IONIZATION OF FOOD AND NUTRACEUTICALS

Robert B. Cody

JEOL USA, Inc., 11 Dearborn Rd. - Peabody, MA 01960 USA

FOr.05-2 **ANCIENT APPLE OF FRIULI VENEZIA GIULIA REGION AS POTENTIAL NEW NUTRACEUTICAL INGREDIENTS**

Stefania Sut⁽¹⁾ - **Stefano Dall'Acqua**⁽²⁾

⁽¹⁾ University of Padova, DAFNAE, Padova - ⁽²⁾ University of Padova, Department of Pharmaceutical and Pharmacological Science, Padova

Keywords: Apple, polyphenols, triterpene, LC-MS, nutraceutical

Introduction

Ancient Apple varieties can present different composition compared to commercial ones [1-6]. For this reason they can be valuable starting material for the development of new nutraceutical ingredients. The aim of this study was to evaluate vitamin C and secondary metabolite contents of nine Ancient Apple of Friuli Venezia Giulia region and compare with commercially available fruits for the possible development of new Apple-based nutraceutical products

Methods

Peels and pulps were separated, dried and extracted in aqueous methanol in ultrasound bath. Hydroxycinnamic acids, catechin, procyanidins, flavonoids, dihydrochalcone and triterpene derivatives composition were determined using LC-DAD-ESI-MSⁿ, LC-ESI-QTOF and LC-APCI-MSⁿ. A targeted metabolomic approach was used to describe the different composition of the ancient varieties compared to commercial ones.

Results

Ancient apple variety showed significant differences compared to commercial varieties with high levels of vitamin C in pulps. Furthermore pulps and peels present higher amount of procyanidins,

catechins and hydroxycinnamic acids. Different pattern in triterpene contents was also observed showing ancient varieties presenting the annurcoic acid, while ursolic and oleanolic characterize the commercial ones. With targeted MS data PLS-DA model was generated and the ancient varieties of FVG region appear to be clusterized mainly on the basis of triterpene and phenolic contents.

Conclusions

Ancient apple can be used for nutraceutical production due to secondary metabolite contents. LC-MS based metabolomics approaches may be useful for description and valorization of these fruits composition. Further studies are needed to optimize extract production from these fruits. The use of ancient apple in nutraceuticals can be useful to protect these fruits threatened by intensive cultivation, underlying the richness of Italian biodiversity.

Novel Aspect

This work point out the usefulness of LC-MS for the characterization of food ingredients for nutraceutical applications and show the valorization of apple biodiversity.

References

1. Boyer, J., Liu, R.H. Apple phytochemicals and their health benefits (2004) *Nutrition Journal*, 3, art. no. 1, pp. 1-45.
2. Belviso, S., Scursatone, B., Re, G., Zeppa, G. Novel data on the polyphenol composition of italian ancient apple cultivars (2013) *International Journal of Food Properties*, 16 (7), pp. 1507-1515.
3. Jakobek, L., Barron, A.R. Ancient apple varieties from Croatia as a source of bioactive polyphenolic compounds (2016) *Journal of Food Composition and Analysis*, 45, pp. 9-15.
4. Donno, D., Galizia, D., Cerutti, A.K. Fruit nutraceutical value in ancient apple cultivars grown in piedmont (Northern Italy) (2012) *ActaHorticulturae*, 940, pp. 131-138.
5. Wojdyło, A., Oszmiański, J., Laskowski, P. Polyphenolic compounds and antioxidant activity of new and old apple varieties (2008) *Journal of Agricultural and Food Chemistry*, 56 (15), pp. 6520-6530.
6. Oszmiański, J., Lachowicz, S., Gławdel, E., Cebulak, T., Ochmian, I. Determination of phytochemical composition and antioxidant capacity of 22 old apple cultivars grown in Poland (2018) *European Food Research and Technology*, 244 (4), pp. 647-662.

FOr.05-3 **RELIABLE IDENTIFICATION OF BIOACTIVE COMPOUNDS IN FOODSTUFFS THROUGH THE USE OF A NOVEL DUAL-FILTER MASS SPECTRAL LIBRARY IN LIQUID CHROMATOGRAPHY**

Francesca Rigano⁽¹⁾ - ***Marianna Oteri***⁽²⁾ - ***Paola Dugo***⁽²⁾ - ***Luigi Mondello***⁽²⁾

⁽¹⁾ *Chromaleont Srl, c/o Dipartimento ChiBioFarAm, University of Messina, Messina* – ⁽²⁾ *University of Messina, Dipartimento ChiBioFarAm, Messina*

Keywords: Linear retention index, bioactive compounds, mass spectral libraries, dual-filter identification process, food samples

Introduction

In the present work, a novel identification approach was developed for the determination of bioactives in foodstuffs by Liquid Chromatography (LC) coupled to mass spectrometry (MS). Non-polar lipids were the analytes of interest, mainly due to their regular LC profile under reversed phase conditions.

Marine organisms centralized special attention because of their high content of ω -3 fatty acids, whose nutraceutical activity was widely studied in the last decades.

Methods

A dual-filter MS/LRI library was built by analyzing more than 30 real-world samples, including vegetable oils, milk and dairy products, sea organisms and food supplements. An odd carbon chain number triacylglycerol series from triononain to trionadecanoin was chosen as the basis of the LRI scale, being these compounds rarely present in real samples and an ultra high performance LC (UHPLC) method, able to maximize lipid separation, was developed.

Results

The MS/LRI contains more than 200 lipid species, including mono-, di- and triglycerides (MG, DG and TG), as well as sterols and sterol esters.

An LRI variability of ± 15 units was obtained for all the compounds, also by changing experimental conditions (column dimension, flow rate, gradient steepness, instrumental setup, temperature), demonstrating the usability of the LRI database at both intra- and interlaboratory levels.

The MS spectrum was essential to discriminate between compounds falling within the same LRI window (completely or partially co-eluted). Under the employed MS conditions, a highly informative fragmentation pattern was obtained. For example, for a TG species, MG and DG fragments were highlighted, so that each MG fragment correspond only to one fatty acid and each DG fragment is related to their combination. The molecule related ion represents the final confirmation of the right fatty acid combination in the TG species.

Conclusions

The novel dual-filter library represents a powerful tool to achieve an univocal identification of non-polar lipids. Further efforts will be required to extend it to other lipid species, such as fatty acids and phospholipids, which represent important components in the clinical application field. Advancements in software, able to perform an automatic identification on the basis of both LRI and MS data, was essential and has been successfully implementing.

Novel Aspect

The complementarity between MS and LRI data was exploited for the first time in LC, allowing to obtain for an unknown sample a very narrow list of candidates for each LC-MS peak.

FOr.05-4 ANTIOXIDANT ACTIVITIES AND PHYTOCHEMICAL PROFILE OF AZORELLA GLABRA WEDD

Immacolata Faraone⁽¹⁾ - *Dilip K. Ray*⁽²⁾ - *Lucia Chiummiento*⁽¹⁾ - *Alka Choudhary*⁽²⁾ - *Daniela Russo*⁽¹⁾ - *Chiara Sinisgalli*⁽¹⁾ - *Eloy Fernandez*⁽³⁾ - *Luigi Milella*⁽¹⁾

⁽¹⁾ Basilicata University, Departement of Science, Potenza - ⁽²⁾ Teagasc Food Research Centre, Department of Food BioSciences, Dublin - ⁽³⁾ Czech University of Life Sciences Prague, Faculty of Tropical AgriSciences, Prague

Keywords: Antioxidant, polyphenols, nutraceutical, Azorellaglabra, Apiaceae

Introduction

The genre *Azorella* belongs to Apiaceae family. The Apiaceae includes different plants commonly used as infusion in South America in folk medicine to treat various diseases such as bronchitis, cold, rheumatism and neuralgia [1]. The aim of this study is to investigate the potential antioxidant properties and the phytochemical profile of the aerial parts of *Azorellaglabra* (*A. glabra*) Wedd, an endemic uninvestigated plant from Bolivia.

Methods

The aerial parts of *A. glabra* were subjected to maceration, solid-liquid, and then to liquid/liquid extraction using different solvents. The polyphenol content, radical-scavenging activity by using

DPPH, ABTS, $\cdot\text{NO}$ and $\text{O}_2^{\cdot-}$ radicals, FRAP and inhibition of lipid peroxidation were determined; and these data were used to calculate Relative Antioxidant Capacity Index (RACI) [2, 3]. The identification and quantification of compounds were performed by LC-ESI-MS/MS.

Results

The ethanol, n-hexane, chloroform, ethyl acetate, n-butanol and water samples were analyzed for their antioxidant activity by different assays. The ethyl acetate fraction reported the highest RACI value (1.17). For this reason, the LC-ESI-MS/MS analysis was focused on the ethyl acetate fraction. The tentative identification and the quantification of compounds present in the ethyl acetate fraction were reached through accurate mass measurements, MS/MS fragmentation patterns and retention times of standards. This fraction revealed the presence of 15 polyphenols and the pentacyclic terpenoid oleanolic acid identified for the first time in *A. glabra*. The most abundant present polyphenols are the orientin and 3,5-di-O-caffeoylquinic acid. These health-promoting compounds have medical and nutraceutical applications [4, 5].

Conclusions

Ethyl acetate fractions of *Azorellaglabra* extract possessed the best antioxidant capacity. The polyphenolic compounds, in particular the abundant orientin and 3,5-di-O-caffeoylquinic acid, may explain the health benefits of the infusions used in folk medicine.

Novel Aspect

This research has evidenced for the first time that *A. glabra* can represent an important natural source for the pharmaceutical, cosmeceutical and nutraceutical applications.

References

1. Tůmová L., Dučaiová Z., Cheel J., Vokřál I., Sepúlveda B., Vokurková D., *Pharmacognosy magazine*, 13(50), 260-264 (2017).
2. Armentano M.F., Bisaccia F., Miglionico R., Russo D., Nolfi N., Carmosino M., Andrade P.B., Valentão P., Diop M.S., Milella L., *BioMedResearch International*, 2015, 11 pages (2015).
3. Dekdouk N., Malafronte N., Russo D., Faraone I., De Tommasi N., Ameddah S., Severino L., Milella L., *Evidence-Based Complementary and Alternative Medicine*, 2015, 9 pages (2015).
4. Praveena R., Sadasivam K., Deepha V., Sivakumar R., *Journal of Molecular Structure*, 1061(5), 114-123 (2014).
5. Beckera L., Zaiter A., Zimmer J.P.D., Karam M.C., Baudelaire E., Scher J., Dicko A., *Industrial Crops and Products*, 87, 116-123 (2016).

POSTER COMMUNICATIONS

MPS-S01 - Nucleic acids

Chairs: Daniele Fabris, Giuliana Bianco

MP-1 / CHARACTERIZATION OF LILBID-MS USING DNA STRUCTURES

Genia Hense - Phoebe Young - Nina Morgner

AK Morgner, Physical and Theoretical Chemistry at the Goethe University Frankfurt, Frankfurt am Main

Keywords: DNA dissociation, LILBID-MS, soft ionization, native MS, non-covalent interaction

Introduction

LILBID-MS is a soft ionization method for mass spectrometry in which non-covalent complexes are transferred into the gas phase via microdroplets of a native-like aqueous buffer, which are irradiated by a mid-IR laser pulse. (1) This leads to an explosive expansion of the sample droplets and release of ions into vacuum. As a result of the desorption process, the ions receive an energy input from the laser pulse, which can be controlled via laser energy transfer parameters.

Methods

This technique is able to detect intact ssDNA and dsDNA. Here we use DNA strands, which consist of 15-20 bp. Under soft laser conditions, we observe nearly no dissociation of the duplex DNA, while harsh conditions lead to the complete dissociation of the dsDNA. We can create melting curves of dsDNA with LILBID by use of a heatable droplet generator. Comparison of LILBID melting curves with UV/Vis shows how LILBID can monitor the solution behavior of DNA.

Results

The tolerance of LILBID-MS towards different buffers and high salt concentrations enables the use of the temperature-independent MgHPO_4 buffer. These buffer conditions ensure the stability of double helices without compromising the spectral resolution. By increasing the temperature of a sample we can achieve different dissociation states of dsDNA. To calculate the duplex ratio, the total peak areas of duplex signals were divided by the sum of all peaks. Plotting the calculated duplex ratio against the temperature of the solution allows the determination of an apparent melting temperature (T_m) of dsDNA because the shape of curves recorded by LILBID-MS mimics the sigmoidal fit of melting curves recorded by UV/Vis spectroscopy. The detected T_m reflects the length of the dsDNAs. Comparison of the apparent LILBID T_m with the measured UV/Vis T_m allows assessing the amount of dissociation caused by energy transfer during the desorption process.

Conclusions

Double-stranded DNAs as non-covalent complexes are suitable for investigating dissociation caused by energy transfer to the molecules during LILBID desorption. Apparent melting curves from LILBID-MS correlate with melting curves recorded by UV/Vis spectroscopy, indicating that LILBID-MS can give reliable information on binding affinities of bimolecular compounds.

Novel Aspect

The method of creating melting curves of dsDNA via LILBID-MS allows the assessment of the amount of

transferred energy during the desorption process.

References

N. Morgner, H. D. Barth, T. L. Schmidt, A. Heckel, U. Scheffer, M. Göbel, P. Fucini, B. Brutschy, *Phys. Chem.* 221, 689-704

MP-2 / THE SIMULTANEOUS DETERMINATION OF INTRACELLULAR NUCLEOSIDE TRIPHOSPHATES AND CYCLIC DINUCLEOTIDES WITH LC/MS/MS

Reinhard Oertel⁽¹⁾ - ***Rayk Behrendt***⁽²⁾ - ***Ali El-Armouche***⁽³⁾

Institute of Clinical Pharmacology, GWT TU Dresden, Dresden (1) - Institute for Immunology, TU Dresden, Dresden (2) - Centrum for Pharmacology and Toxicology, TU Dresden, Dresden (3)

Keywords: dNTP/NTP, dinucleotide, LC/MS/MS, isomeric separation

Introduction

The determination of nucleoside triphosphates and cyclic dinucleotides in cells is of interest for various issues: The cyclic dinucleotides adenosine monophosphate-guanosine monophosphate (cGAMP) and diGMP could be a tumor marker in human fibroblasts. The quantification of deoxynucleoside triphosphates (dNTPs) and the corresponding ribonucleotide triphosphates (NTPs) can be used to measure enzyme activity in knockout cells

Methods

ATP and dGTP are isomers, CTP, UTP and dTTP have very similar masses. Therefore a complete chromatographic separation was an aim of the study. An HPLC-Method with a reversed phase column and a mobile phase gradient were applied to separate all analytes in less than 10 min. Two transitions per substance were used in multiple reaction monitor (MRM) mode of the tandem mass spectrometer.

Results

An 8-point external calibration curve in buffer at the range from 7.8 – 1000 ng/mL for 10 nucleoside triphosphates and cyclic dinucleotides was prepared and measured with samples of the assay. Quantification was performed by peak area method. A weighted (1/x) regression first order was performed to determine the concentration of the analytes. The calibration curves were accurate and precise over the range from 7.8 – 1000 ng/mL and the correlation coefficient (r) exceeded 0.996 in all cases. The concentrations 7.8 and 1000 ng/mL are the lower and the upper limit of quantification (LLOQ and ULOQ). Coefficients of variation and accuracy were calculated for each standard and are supposed to vary by less than 15%. All tested real samples contained nucleotide. To estimate the matrix effect, the samples were additionally analyzed by the standard addition method.

Conclusions

A simple and sensitive LC / MS / MS method was developed to separate different, partially isomeric nucleotides. The method provides a convenient tool for measuring intracellular dNTP / NTP levels that can be used to measure the molecular phenotype of mouse cells. In addition, cyclic dinucleotide concentrations in human cells were determined to be potential tumor markers.

Novel Aspect

The novel aspect of the study is that with only one measuring method different nucleotides in different matrices can be determined simultaneously for different questions.

MP-3 / SEMI-QUANTITATIVE DETERMINATION OF OLIGONUCLEOTIDE DRUG IMPURITIES: MAIN AND CO-ELUTING SPECIES

Stilianos G. Roussis - Claus Rentel

Ionis Pharmaceuticals, ADQC, Carlsbad

Introduction

An advanced Ion Pair Reverse Phase (IP RP) LC MS characterization method for impurities co-eluting with the main component is currently available and extensively used in our labs for release and stability testing of oligonucleotide therapeutics. However, the method requires individual extraction of each impurity with integration of the resulting ion chromatogram. For certain applications, *e.g.* optimization of synthesis conditions, purification development or evaluation of reagents, the number of samples requiring analysis is often large, but a semi-quantitative determination comparing impurity levels between the samples is sufficient. In this work, a process is described for using the average mass spectrum of the main product LCMS peak to automatically determine the amounts (%area) of the co-eluting impurities.

Methods

Ion Pair-Reversed Phase (IP-RP) HPLC chromatography coupled to low- and high-resolving power mass spectrometers has been used for the experiments using Agilent instrumentation. For an oligonucleotide drug component a list of chemical formulas of known impurities is generated using an established template. An algorithm was developed to automatically determine the isotopic distribution of each formula and its most abundant mass (MAbM). For each chemical formula, a list of MAbMs is determined corresponding to the charge states of the impurity in the ESI spectrum. An average mass spectrum is obtained by summation of the signal across the main LC peak incorporating the trace impurities. The relative amounts of the co-eluting impurities are determined by integration of the peak abundances of each species determined in the spectra.

Preliminary Data

The method has been extensively applied for the automated comparison of co-eluting oligonucleotide impurity levels in samples. For five different oligonucleotides, results have also been compared to those obtained by the quantitative conventional method, and were in very good agreement. The average absolute difference between the two approaches was 0.3%, with the best agreements achieved when the square root of the measured ESI abundance was used for the calculations. The repeatability of the method was evaluated, indicating an average standard deviation of 0.04% (calculated from all detected species). Application of the method for monitoring impurity changes in different lots of drug product readily identified the specific impurities statistically deviating from established values. Principal component analysis (PCA) was used to determine outliers. The loadings of PCA analysis were used to classify a set of 10 oligonucleotides in classification groups according to similarities in the formation of the same type of impurities. A plot of the PCA scores outlined the impurity types primarily contributing to the classification. For some impurities, multiple linear regression revealed a linear dependence of the total amounts detected to the number and type of a combination of base residues in the oligo. The current approach due to its rapidity and simplicity can be readily applied to determine relationships between the impurities detected and the experimental process development and manufacturing conditions potentially affecting them.

Novel Aspect

Development and application of an automated approach for the determination of impurity species co-eluting with the main oligonucleotide drug component.

MP-4 / INTERACTIONS OF OXALIPLATIN AND CARBOPLATIN WITH DNA BUILDING BLOCKS AS STUDIED IN THE GAS PHASE: A COMBINED EXPERIMENTAL AND THEORETICAL STUDY

Jean-Yves Salpin⁽¹⁾ - *Maria Elisa Crestoni*⁽²⁾ - *Davide Corinti*⁽²⁾ - *Simonetta Fornarini*⁽²⁾ - *Debora Scuderi*⁽³⁾ - *Barbara Chiavarino*⁽²⁾

CNRS, LAMBE - UMR 8587 - Université d'Evry Val d'Essonne, Evry (1) - Università di Roma la Sapienza, Dipartimento di Chimica e Tecnologia del Farmaco, Roma (2) - Laboratoire de Chimie-Physique (LCP), Université Paris-Sud, Orsay (3)

Keywords: Pt-based anti-tumor agents, IRMPD spectroscopy, DFT calculations, nucleic acid building blocks

Introduction

In spite clinical successes for treating different types of tumors, cisplatin [cis-Pt(NH₃)₂(Cl₂)] displays severe toxic side effects and therapy resistance. This has motivated searches for structurally and/or functionally analogous alternatives, leading to carboplatin and oxaliplatin as second and third-generation platinum drug. The aim of this work is to study the interactions occurring in the gas phase between these two drugs and a set of DNA building blocks.

Methods

10⁻⁴ M water solutions of either carboplatin or oxaliplatin and analyte (5'-dGMP) or (5'-dAMP) were mixed together and stored at room temperature for 24 h, to allow the formation of the different complexes. Complexes were generated in the gas phase by electrospray ionization. Their structure was probed by combining mass spectrometry and IRMPD spectroscopy to DFT calculations at the B3LYP/6-311G(d,p) level.

Results

This combination of approaches allows operating in the gas phase, which makes it possible to overcome any solvent or counter-ion effect, and thus allows describing the intrinsic reactivity of the metal cations. Presently, the complexes generated in the gas phase with adenine(A) and guanine (G) nucleobases are of the type [(platinum drug)H(nucleobase)]⁺. Our DFT survey indicates that protonation preferentially occurs on one of the carbonyl group of both oxaliplatin and carboplatin. This is consistent with mechanisms proposed in previous studies. [1,2] Comparison between IRMPD spectra and DFT calculations suggests that both oxaliplatin and carboplatin specifically interact with the N7 position of guanine. This interaction mode is therefore similar to what was observed with cisplatin. [3] Data obtained with adenine are more delicate to interpret. IRMPD spectroscopy points to a mixture of structures of comparable relative energies with both platinated agents, corresponding to different coordination schemes.

Conclusions

The present study shows that the interactions occurring between these two platinum drugs are similar to those observed with cisplatin. Concerning guanine, the interaction specifically occurs at the N7 position, in agreement with the remarkable stability of the corresponding structures. As for adenine, interaction with either N1, N3 or N7 sites results in structures which are close in energy, in agreement with the mixture of forms observed experimentally.

Novel Aspect

To our knowledge, this work is the first gas-phase experimental study of the interactions taking place between carboplatin and oxaliplatin, with A and G.

References

1. Lucas M. F. A., Pavelka M., Alberto M. E., Russo N., J. Phys. Chem. B, 113, 831 (2009).
2. Alberto M. E., Lucas M. F. A., Pavelka M., Russo N., J. Phys. Chem. B, 113, 14473 (2009).
3. Chiavarino B., Crestoni M. -E., Fornarini S., Scuderi D., Salpin J-Y., J. Am. Chem. Soc., 135, 1445 (2013).

MP-5 / DEVELOPMENT OF THE LC-MS METHOD FOR THE ANALYSIS AND CHARACTERIZATION CHEMICALLY MODIFIED RNA

Dominika Strzelecka⁽¹⁾ - *Miroslaw Smietanski*⁽²⁾ - *Marcin Warminski*⁽¹⁾ - *Pawel Sikorski*⁽²⁾ - *Jacek Jemielity*⁽²⁾ - *Joanna Kowalska*⁽¹⁾

Faculty of Physics, University of Warsaw, Warsaw (1) - Centre of New Technologies, University of Warsaw, Warsaw (2)

Keywords: RNA, thiophosphate, modified nucleotide, isotope labeling, quantitative analysis

Introduction

Messenger RNA (mRNA) is a promising therapeutic agent with potential uses in cancer therapy, gene therapy, and cell reprogramming. The stability of mRNA largely depends on the 3'- end and 5'- end structures, which are targets for enzymatic degradation. Chemical modifications in those regions can be introduced into in vitro transcribed mRNAs in order to confer on them resistance to degradation and increase their translation efficiency, thereby increasing their therapeutic potential.

Methods

We developed an LC-MS method for quantitative assessment of incorporation of ATP analogs modified at the α -phosphate (ATP α S, ATP α BH₃) into transcripts by RNA polymerases. Modified RNA was enzymatically degraded to single nucleotides, which were then resolved and quantified by ion pair chromatography coupled with ESI-QQQ. In order to improve the quantification, we synthesized RNA degradation products labelled with heavy oxygen (¹⁸O) within the phosphate and used them as internal standards.

Results

We obtained chemically modified RNA by in vitro transcription, wherein standard ATP was replaced or mixed at different ratios with an appropriate chemically modified ATP analog (ATP α X) containing an O-to-X substitution (X=S, BH₃) at the alpha phosphate and differing in stereochemistry (SP versus RP). The assessment of efficiency and frequency of incorporation of the modified NTP versus natural NTP by SP6 and T7 RNA polymerase was analyzed by LC-MS. First, RNA was subjected to complete degradation into single 5'-nucleotides by SVPDE enzyme. The concentrations of AMP and phosphate-modified AMP analog (AMPX) in the resulting mixture was quantified by LC-MS method employing ¹⁸O-labelled standards. The analysis of AMP:AMPX ratio in the analyzed transcript as a function of ATP:ATP α S ratio in the transcription reaction provided insights into the structure of RNA as well as substrate stereoselectivity of RNA polymerases.

Conclusions

Proposed LC-MS quantitative method enabled determination of frequency of incorporation of different modified nucleotides into in vitro transcribed RNA, thereby it providing a tool for establishing mRNA structure-biological properties relationship.

Novel Aspect

The proposed procedure could contribute to popularization of the application of low resolution equipment in the analysis of large molecules such as nucleic acids.

Acknowledgement: This work is supported by the National Science Centre, Preludium (UMO-2016/23/N/ST4/03186)

MP-6 / SOLUTION AND GAS-PHASE STRUCTURES OF DNA G-QUADRUPLEXES: ION MOBILITY MASS SPECTROMETRY, SOLUTION SPECTROSCOPY, AND COLLISION CROSS SECTION CALCULATIONS

Anirban Ghosh - *Frédéric Rosu* - *Valérie Gabelica*

INSERM, University of Bordeaux/IECB, Pessac, France

Keywords: G-quadruplex, conformational polymorphism, Ion mobility mass spectrometry, Collision Cross section, $^1\text{H NMR}$

Introduction

G-quadruplexes are guanine-rich sequences mostly found in promoters and telomeres associated with key cellular processes [1-2]. They are polymorphic in terms of number of tetrads, strands orientation, cation involved and loop arrangement[3]. Here we are using Ion mobility mass spectrometry (IM-MS) to characterize conformations of diverse sequence topology in native MS buffer [4]. This is further supplemented by $^1\text{H NMR}$, CD and UV melting.

Methods

All IM-MS experiments were performed on Agilent 6560 IMS-Q-TOF in negative ion mode with Helium in the drift tube [5]. All circular dichroism (CD) and UV melting experiments were carried out on Jasco J-815 and SAFAS UV mc² spectrophotometer at 25°C. All $^1\text{H NMR}$ experiments were performed on Bruker Avance 700 MHz with TXI probe. The theoretical CCS values for each charge states were calculated from Gaussian optimized conformations by modified EHSS method by Siu et al. [6-7].

Results

CD in native MSbuffer (100 mM TMAA + 1 mM KCl) shows preservation of the high-[KCl] G-quadruplex conformation. This is further supported by the $^1\text{H NMR}$ spectra of the imino proton (10-14 ppm) regions, which agree well with the literature. ESI-MS shows a major population of 2-K⁺ complex for all 3-quartet conformations, and also a 1-K⁺ complex (2-quartet) for a few human telomere sequences. The collision cross-section (CCS) distributions are broad for all hybrid and anti-parallel conformations, and narrower for parallel conformations. For a given sequence length and charge state, the loop configuration, quartet polarity, bulges and flanking bases give rise to significant differences in CCS. The starting PDB structures give a higher theoretical CCS values compared to experimental ones. The structure neutralized to match the experimental charge states are optimized in the gas phase at the semi-empirical PM7 level and closely match the experimental CCS values. EHSSrot (Siu parameters) or the TM are used to calculate the theoretical CCS values.

Conclusions

Most examined G-quadruplexes fold with the same topology in native MS (1 mM KCl) and traditional high KCl content buffers. IM coupled with ESI-MS is able to distinguish different shapes formed in solution. With these well-characterized G-quadruplexes, we validate the gas-phase structure optimization with PM7 to generate 3D models from which to calculate the CCS values

Novel Aspect

G-quadruplexes fully characterized by CD, UV and NMR in MS-compatible solutions were used to validate the interpretation of IM-MS data in terms of structure.

References

1. Maizels, N. Gray, L. T., PLoS Genet., 9, e1003468 (2013).
2. Rhodes, D., Lipps, H. J., NAR, 43, 8627-8637 (2015).
3. Dai, J., Carver, M., Yang, D., Biochimie, 90, 1172-1183 (2008).
4. Marchand, A., Gabelica, V., JASMS, 25, 1146-1154 (2014).
5. Marchand, A., Livet, S., Rosu, F., Gabelica, V., Anal. Chemistry, 89, 12674-12681 (2017).
6. Shvartsburg, A. A., Mashkevich, S. V., Baker, E. S., Smith, R. D., JPC A, 111, 2002-2010 (2007).
7. Siu, C. -K., Guo, Y., Saminathan, I. S., Hopkinson, A. C., Siu, K. M., JPC B, 114, 1204-1212 (2009).

MP-7 / CLUSTER-ASSEMBLED SMART SURFACES FOR ON-TARGET PROCESSING OF BIOLOGICAL SAMPLES IN MALDI-MS ANALYSIS

Emanuele Barborini⁽¹⁾ - ***Simone Vinati***⁽¹⁾ - ***Giacomo Bertolini***⁽¹⁾ - ***Marc Baumann***⁽²⁾

Tethis SpA, R&D, Milan (1) - University of Helsinki, Clinical Proteomics Unit, Helsinki (2)

Keywords: MALDI target, Imaging, processing, film, nanoparticles

Introduction

Matrix-assisted laser desorption/ionization mass spectrometry (MALDI-MS) is a powerful analytical tool endowed with the capability to analyze biological samples ranging from biomolecules to micro-organisms, up to whole histological sections. Due to the complexity of such analytes, there is the need for ad-hoc sample pre-processing (*e.g.* washing, reduction/alkylation, enzymatic cleavage), which can be performed on-target, through the use of smart surfaces.

Methods

Cluster Beam Deposition [1,2] was adopted to deposit patterned functional films made of ultrafine TiO₂ nanoparticles, either on polymeric, metallic and glass slides at room temperature. Peculiar nanoporosity and bio-affinity of this material play a synergic role in sample capturing. TiO₂ super-hydrophilicity by UV irradiation [3] provides uniform spreading of liquid samples and optimal adhesion of histological sections onto the patterned functional areas.

Results

Hydrophobic barrier at the border of the super-hydrophilic functional areas acts as an effective droplet confinement structure, which allows easy and reliable on-target processing of liquid-phase samples. In comparison with standard in-vial approach, on-target processing of protein samples (including secondary cleaning and Trypsin digestion) avoids the loss of sample fractions, due to cleaning, as well as of peptide fractions due to their capture on vial plastic walls. Regarding histological samples, improved adhesion avoids detachment, tearing and shrinking, during tissue processing (*e.g.* dehydration, delipidation, fixation). Improved adhesion can also benefit pharmaco-kinetic studies, where any tissue treatment should a-priori be excluded [4]. Tissue stability permits matrix removal after MALDI Imaging, in order to carry out immunohistochemical analysis on the very same section. In general, nanoscale surface roughness offers distributed nucleation points favoring uniform and homogeneous crystallization of MALDI matrix.

Conclusions

Advanced surface engineering concepts have been adopted to develop innovative MALDI targets aiming to reliable on-target sample processing. Advantages were thoroughly demonstrated with protein samples, subjected to secondary cleaning and Trypsin digestion, and with histological samples subjected to harsh chemical treatments prior MALDI Imaging. Benefits are also expected in pathogen analysis as well as in robotic sample processing for clinical diagnostics.

Novel Aspect

Improved on-target sample processing in MALDI-MS through the use of nanoporous bio-affine films made of ultrafine TiO₂ nanoparticles deposited on any kind of a substrate.

References

1. Wegner K. et al. Journal of Physics D: Applied Physics, 39 (2006) R439.

2. Barborini E. et al. European Physical Journal D, 24 (2003) 277.
3. Wang R. et al. Nature 388 (1997) 431.
4. Calza G. et al. Frontiers in Pharmacology, 9:141 (2018) 1.

MP-8 / PROTEOME PROFILING OF BRONCHIOALVEOLAR LAVAGE FLUID AFTER EXPOSURE TO MAGNETIC IRON AND COBALT OXIDE NANOPARTICLES

Anja M. Billing⁽¹⁾ - **Kristina B. Knudsen**⁽²⁾ - **Ulla Vogel**⁽³⁾ - **Frank Kjeldsen**⁽¹⁾

University of Southern Denmark, Department of Biochemistry and Molecular Biology, Odense (1) - National Research Centre for the Working Environment, National Research Centre for the Working Environment, Copenhagen (2) - Technical University of Denmark, Department of Micro- and Nanotechnology, National Research Centre for the Working Environment, Copenhagen (3)

Keywords: clinical proteomics, label-free, EvoSep, bronchioalveolar lavage, iron cobalt oxide nanoparticles

Introduction

Magnetic nanoparticles are of great interest in biomedicine (e.g. magnetic resonance imaging, targeted tumor therapy), however, their biocompatibility is an ongoing concern and large-scale in vivo toxicology studies are mostly lacking. Here we examined the effects of magnetic iron and cobalt oxide nanoparticles on the lung in mice using label-free quantitative proteomics. The new EvoSep LC system, specifically developed for clinical proteomics, was evaluated.

Methods

Toxicity of iron cobalt oxide nanoparticles was examined in vivo by intratracheal instillation in mice. Bronchioalveolar lavage fluid (BALF) samples were collected at day 1 and 3 post exposure. Nanoparticles with different iron cobalt composition were included (Fe Co: 1:0, 3:1, 1:1, 1:3, 0:1) testing two doses (54 µg, 162 µg). The BALF proteome of 168 mice was profiled by label-free LC-MS/MS (EvoSep coupled to Orbitrap Fusion Lumos, Thermo) using 21 min gradients.

Results

Performance of EvoSep was evaluated. Average Pearson correlation coefficient of 0.99 was determined among 16 technical replicates demonstrating very high reproducibility. To increase identification and quantification a mix of representative samples was separated by high pH fractionation and measured in parallel. Biological replicates (168 BALF samples) were run as single shot injections using 21 min gradients. A total of 1400 proteins were quantified with average of 400 proteins per biological replicate (min: 269, max: 602). Statistical analysis was performed by empirical Bayes moderate T test incorporated in the Limma package. Differentially expressed proteins were clustered and enrichment analysis was performed. Highly enriched were proteins associated with extracellular vesicles, stress and immune response, as well as glutathione metabolism, response to metal ion and reactive oxygen species. To assist pathway identification in-depth BALF library from high pH fractionation was included in the bioinformatics analysis.

Conclusions

We present here a proteomic toxicity screening testing simultaneously 6 different nanoparticle types (magnetic iron cobalt oxide nanoparticles with Fe Co ratio: 1:0, 3:1, 1:1, 1:3, 0:1; carbon black as control), after pulmonary inhalation at two time points and two concentrations. As demonstrated EvoSep combined with short LC gradients and high-resolution mass spectrometry enables high-throughput proteome profiling of clinical samples such as BALF or plasma.

Novel Aspect

The EvoSep LC system designed for clinical proteomics allowing for high throughput and robust

measurement was tested on real-life samples.

MP-9 / CONCURRENT QUANTIFICATION OF SYNAPTOTAGMIN-1 AND SNAP-25 IN CLINICALLY RELEVANT VOLUMES OF CSF BY IMMUNOAFFINITY LC-MS/MS

Ann Brinkmalm⁽¹⁾ - *Annika Öhrfelt*⁽¹⁾ - *Lutz Frölich*⁽²⁾ - *Lucrezia Hausner*⁽²⁾ - *Henrik Zetterberg*⁽¹⁾ - *Kaj Blennow*⁽¹⁾

The Sahlgrenska Academy at University of Gothenburg, Institute of Neuroscience and Physiology, Department of Psychiatry and Neurochemistry, Mölndal (1) - University of Heidelberg, Department of Geriatric Psychiatry, Central Institute for Mental Health Mannheim, Mannheim (2)

Keywords: PRM-MS, immunoprecipitation, cerebrospinal fluid, synaptic proteins

Introduction

Reduction of synaptic activity is an early event in several neurodegenerative diseases. In two earlier studies we have used affinity purification combined with PRM-MS to demonstrate increased CSF levels of the synaptic proteins synaptotagmin-1 and SNAP-25 in patients with Alzheimer's disease compared with controls(1, 2). In this study the strategy was further developed to concurrently measure levels of both proteins in individual patient CSF samples.

Methods

The immunoprecipitation was performed on a KingFisher™ Flex System which uses magnetic rods to move particles through the various binding, mixing, washing and elution phases in a 96 well plate format. HR-PRM analyses were performed on a Q Exactive™ Hybrid Quadrupole-Orbitrap™ MS coupled to an Ultimate 3000 standard liquid chromatography system. The instrument was set to acquire scheduled pairs of PRM scans in profile mode.

Results

In this study, we developed and validated a high-throughput assay for concurrent quantification of synaptotagmin-1 and SNAP-25 in clinically relevant volumes of CSF. The assay was utilized to compare the levels of synaptotagmin-1 and SNAP-25 in a clinical cohort with CSF samples from patients diagnosed with early AD, prodromal AD, mild cognitive impairment, late AD, mixed AD, and non-demented controls. The increased CSF levels of both synaptotagmin-1 and SNAP-25 in AD patients compared to controls are in agreement with the results from two previous studies where CSF synaptotagmin-1 and SNAP-25 were quantified separately.

Conclusions

We have shown that automated immunoprecipitation combined with PRM-MS is robust and sensitive enough to concurrently measure two synaptic proteins in CSF from individual patients.

Novel Aspect

This platform will make it feasible to perform larger and diagnostic diversified studies as well as adding other biomarkers mirroring other pathologies.

References

1. Brinkmalm A, Brinkmalm G, Honer WG, Frölich L, Hausner L, Minthon L, et al. SNAP-25 is a promising novel cerebrospinal fluid biomarker for synapse degeneration in Alzheimer's disease. *Molecular Neurodegeneration*. 2014;9.
2. Öhrfelt A, Brinkmalm A, Dumurgier J, Brinkmalm G, Hansson O, Zetterberg H, et al. The pre-synaptic vesicle protein synaptotagmin is a novel biomarker for Alzheimer's disease. *Alzheimers Res Ther*.

MP-10 / SINGLE SHOT DIA PROFILING OF >1500 PLASMA PROTEOMES OF THE WEIGHT LOSS AND MAINTENANCE STUDY DIOGENES

Roland Bruderer⁽¹⁾ - **Jan Muntel**⁽¹⁾ - **Sebastian Müller**⁽¹⁾ - **Oliver M. Bernhardt**⁽¹⁾ - **Tejas Gandhi**⁽¹⁾ - **Polina Mironova**⁽²⁾ - **Ondine Walter**⁽²⁾ - **Jérôme Carayol**⁽²⁾ - **Arne Astrup**⁽³⁾ - **Wim H. M. Saris**⁽⁴⁾ - **Jörg Hager**⁽²⁾ - **Armand Valsesia**⁽²⁾ - **Loïc Dayon**⁽²⁾ - **Lukas Reiter**⁽¹⁾

Biognosys AG, RnD, Schlieren (1) - Nestle, Nestlé Institute of Health Sciences, Lausanne (2) - University of Copenhagen, Department of Nutrition, Exercise and Sports, Copenhagen (3) - Maastricht University, Maastricht University Medical Centre, Maastricht (4)

Introduction

Comprehensive, robust, high-throughput analysis of the plasma proteome has the potential to enable analysis of the health state. Up to now, analysis of the plasma proteome was achieved using delicate nano-flow setups. To reduce this limitation, we established a robust capillary-flow LC-MS-DIA setup capable of measuring 32 plasma proteomes a day. Using this setup, we acquired the samples of a large-scale weight loss and maintenance plasma sample study (DiOGenes).

Methods

Plasma samples were randomized and prepared with an optimized in-solution digestion protocol and spiked with Biognosys' Plasma Deep Dive. The samples were acquired on a Thermo Scientific Orbitrap Fusion Lumos connected to a Waters M-class UPLC and a 300 µm*150 cm Waters CSH 1.7 µm column. DDA runs were performed using HPRP fractionation. DIA was done in a 45 min injection-to-injection mode. Analysis of DIA runs was performed using Spectronaut Pulsar.

Results

After testing several combinations, we obtained optimized conditions using a 300 µm inner diameter column. This setup has the advantage of offering increased flow rates (5 µl/min) resulting in small injection-to-injection overhead times (*i.e.*, 5 min). A gradient of 40 min was determined as an optimal gradient length, resulting in 90% of the maximally achieved identifications with longer gradients. Further characterization of the system showed precise quantification of 47 FDA approved plasma biomarkers (CVs below 20%).

We successfully used this capillary LS-MS DIA setup to acquire 1,508 (plus 66 control pools) samples of the DiOGenes project. We identified in average 450 proteins per samples and achieved a dataset completeness of 77% at the protein level. Proteins of the inflammation system were the most variable and proteins involved in blood coagulation were the least variable protein groups. Unsupervised clustering showed that the weight loss was most severe intervention accompanied with significant protein changes in plasma.

Conclusions

Differential abundance analysis identified proteins to be significantly differential between all of the time points revealing biology. Non-enzymatic addition of glucose to proteins (*i.e.*, glycation) was detected and could be followed over the time course of the study. The Plasma Deep Dive kit could be used to perform extrapolation to all identified proteins enabled label-free absolute quantification of 490 additional proteins in of 1,472 plasma samples.

Novel Aspect

High-throughput capillary-flow DIA of 1,574 human plasma proteomes with label-free absolute quantification and unbiased determination of glycation sites.

MP-11 / IDENTIFICATION OF COMPLETE POOL OF BIO ACTIVE PEPTIDES IN MOUSE BRAIN STRIATUM INCLUDING SEPS, A NOVEL CLASS OF NON-CLASSICAL BIO ACTIVE PEPTIDES

Harshavardhan Budamgunta

Center for proteomics, University of Antwerp/Biology, Antwerp

Keywords: short Open Reading Frame (sORF), Mass spectrometry, Bio active peptides, Peptidomics, sORF encoded polypeptides (SEPs).

Introduction

Bio active peptides are involved in the regulation of major physiological processes that occur in an organism [1]. Classical bio active peptides are cleaved from their precursor and are directed towards the secretory pathway [2]. Recently, another novel class of bioactive peptides often termed as SEPs (sORF encoding polypeptides) alias micro peptides alias non-classical bio active peptides were discovered [3]. These are called non-classical because they are translated directly from small ORFs [3]. Therefore, they lack a signaling sequence and are released into the cytoplasm, not directed through the secretory pathway [3]. Despite their unconventional expression, some of these peptides are found to play critical functional roles [4-10]. The first SEP was discovered in *C. Elegans* in 1996 [11]. From then, many sORFs were discovered in plants such as Polaris gene (*Arabidopsis thaliana*), ROTUNDIFOLIA4, Brick, Zm401 and Zm908p1 [12-16]. Recently, with the advent of ribosome profiling (an RNA sequencing technique), millions of sORFs potentially translated to peptides were discovered in *Saccharomyces cerevisiae*, *Drosophila melanogaster*, *Mus musculus*, *Danio rerio* and *Homo sapiens* by various research groups [17-18]. Although ribosome profiling experiments predicted SEPs that are being translated, proving their actual bio chemical existence is quite challenging. In 2012, with the aid of mass spectrometry 90 SEPs were identified in a human cell line [19]. This was followed by the identification of 237 SEPs in humans in 2017 [20]. Since their discovery, millions of SEPs were predicted through ribosome profiling data but only few (hundreds) SEPs were biochemically identified. Four challenges hamper these identifications. The peptides are thought to be in Low abundance (1) – In general natural peptides are less abundant than cytosolic proteins; 2) Narrow mass range – The mass range of SEPs is between 0 -10 KDa which is often missed in the regular proteomic analysis, therefore they are absent in most proteomic data repositories. In addition they often lack of cleavage site for the proteases; 3) Difficulties in data analysis – Search engines are typically not good to identify long peptides with multiple charges; 4) False positives – (During the database search) The databases that are typically used are concatenated ones that contains both known large proteins and predicted SEPs. Due to (up to) a 100 fold larger number of the predicted SEPs (millions) in the database when compared to known uniprot entries, false discoveries are likely to aggregate more in the SEP database. In present study, we performed an extensive analysis on mouse brain striatum sample with an intention to identify SEPs. In order to do so, we developed a protocol to overcome the challenges that arise during the process of identifying these peptides with a high confidence.

Methods

The mouse brain peptides were extracted using a mixture of methanol:water:acetic acid (90:9:1). The peptides were fractionated into 20 fractions in reverse phase liquid chromatography. All the fractions are analyzed in nanoLC QExactive plus orbitrap mass spectrometer, 60 min linear gradient for each fraction. Two replicates for each fraction were run.

Peptide identification: A concatenated database was used: “mouse uniprot reviewed proteome downloaded on 2017-03-01 + predicted SEPs from mouse ribosome profiling data (sorfs.org) + contaminant database [2017-03-01]”. The database searching is performed in three methods. PEAKS (version 7.0, Bioinformatics Solutions Inc.) was used as the first strategy for the database search. Variable modifications used in the settings are: oxidation of methionine, C-terminal amidation of peptides, N-terminal acetylation of proteins, pyro-glu formation from glutamine and glutamic acid at N-terminus, acetylation of lysine, phosphorylation of tyrosine, threonine and serine.

The peptides are searched with Mascot (Version 2.1, Matrix science Pvt Ltd) and Sequest HT (University of Washington, USA) with no enzyme selected. The spectra from precursors with mass range from 500 to 10000 Daltons were selected for the identification process. The database search was conducted on deconvoluted data that was exported from PEAKS. The precursor mass tolerance was set to 10 ppm and the fragment ion mass tolerance was set to 0.02 Daltons. Only high confident peptides (1% False Discovery Rate) were chosen in this analysis. Dynamic modifications including N-terminal acetylation of protein, N-terminal pyro-glutamic acid formation from glutamine, C-terminal amidation of peptide and oxidation of methionine were set as variable modifications. A global false discovery rate (FDR) for peptide and protein identification was set to 1%.

The third approach is similar to the second approach but without using a deconvolution step.

Peptide validation: After identifying peptides with high confidence. The resulting peptides were validated through four steps.

Manual verification of the fragmentation spectrum for the ion series match.

BLASTp analysis for the novelty of the identified SEPs. In this section, the fragment that is matched with the in-silico generated ion series was only searched through BLASTp but not the whole identified sequence.

MS2PIP based validation: MS2PIP is a tool that generate predicted ion peak signal intensities based on artificial neural network. The experimentally identified fragmentation spectra was be matched to the MS2PIP generated fragment intensities [21].

The sORF attributes such as ORF score, FLOSS score, Ribosome occupancy, conservation, and in frame coverage of the identified SEPs will be checked [18,22].

Only those peptides passing all validation steps, are considered confident. Finally experimental fragmentation spectra of the identified peptides were compared to that obtained from a synthetic equivalent.

Results

In total 5650 unique peptides were identified combining both the replicates. Using charge state deconvolution of peptides, we observed 36% increase in unique peptide identifications and we identified longer peptides ranging from 0-85 AA long. Classical bio active peptides: 926 peptides originating from 37 neuropeptide precursors were identified. Non classical bio active peptides: 12 peptides originating from 4 SEPs were identified.

Prior to the validation steps, 53 putative SEPs were identified with 1% FDR (this number can be comparable to the previous studies identifying 90 SEPs). Out of the 53, 49 identifications passed through the first filter (manual sequence verification). Out of 49, only 11 identifications passed through the second filter (BLASTp based verification). Out of 11, 6 identifications survived the third filter (MS2PIP based validation). Out of 6, 4 SEPs passed through the fourth filter (sORF attributes cross check).

Conclusions

A protocol to identify a complete set of natural bioactive peptides is created. Using this approach we identified 926 peptides originating from 37 neuropeptide precursors. To our knowledge, using peptidomics approach, this is the highest number of peptides derived from neuropeptide precursors identified in mouse brain in any study so far. Using a four step validation approach we effectively reduce the chances for eliminated false positive identifications and confidently identified 4 SEPs, out of which 3 SEPs are novel.

Novel Aspect

First, identification of 926 peptides originating from 37 neuropeptide precursors through peptidomics approach in mouse brain striatum is the highest number of unique peptide identifications reported till date. Second, the concept of SEPs is relatively novel which originated in 1996. There are studies proving the functional importance of the SEPs. These SEPs can be a heavy interest for the pharmaceutical companies. Most importantly, we observed false positive identifications in the data where few studies claim to identify SEPs. In this study we proposed a method to eliminate the false positives and result with confidently identified SEPs. In our study we confidently identified 4 SEPS, out of which 3 SEPS are novel and not reported previously.

References

1. C. Tamvakopoulos, *Mass Spectrom Rev* 2007, 26, 389.
2. L. D. Fricker, *Mol Biosyst* 2010, 6, 1355.
3. J. Ruiz-Orera, P. Verdaguer-Grau, J. L. Villanueva-Cañas, X. Messeguer, M. M. Albà, *bioRxiv* 2017, 10.1101/064915.
4. P. Yang, C. Read, R. E. Kuc, G. Buonincontri, M. Southwood, R. Torella, P. D. Upton, A. Crosby, S. J. Sawiak, T. A. Carpenter, R. C. Glen, N. W. Morrell, J. J. Maguire, A. P. Davenport, *Circulation* 2017, 135, 1160; S. C. Chng, L. Ho, J. Tian, B. Reversade, *Dev Cell* 2013, 27, 672; A. Pauli, M. L. Norris, E. Valen, G. L. Chew, J. A. Gagnon, S. Zimmerman, A. Mitchell, J. Ma, J. Dubrulle, D. Reyon, S. Q. Tsai, J. K. Joung, A. Saghatelian, A. F. Schier, *Science* 2014, 343, 1248636.
5. C. Lee, K. Yen, P. Cohen, *Trends Endocrinol Metab* 2013, 24, 222.
6. A. Matsumoto, A. Pasut, M. Matsumoto, R. Yamashita, J. Fung, E. Monteleone, A. Saghatelian, K. I. Nakayama, J. G. Clohessy, P. P. Pandolfi, *Nature* 2017, 541, 228.
7. N. G. D'Lima, J. Ma, L. Winkler, Q. Chu, K. H. Loh, E. O. Corpuz, B. A. Budnik, J. Lykke-Andersen, A. Saghatelian, S. A. Slavoff, *Nat Chem Biol* 2017, 13, 174.
8. A. A. Nakamura, R. ; Mukai, M. ; Kobayashi, S. ; Lasko, P. F., *Science* 1996, 276, 2075.
9. J. Zanet, H. Chanut-Delalande, S. Plaza, F. Payre, *Curr Top Dev Biol* 2016, 117, 199; J. Savard, H. Marques-Souza, M. Aranda, D. Tautz, *Cell* 2006, 126, 559.
10. D. M. Anderson, K. M. Anderson, C. L. Chang, C. A. Makarewich, B. R. Nelson, J. R. McAnally, P. Kasaragod, J. M. Shelton, J. Liou, R. Bassel-Duby, E. N. Olson, *Cell* 2015, 160, 595; B. R. M. Nelson, C. A. ; Anderson, D. M. ; Winders, B. R. ; Troupes, C. D., *Science* 2016, 271.
11. K. P. Van de Sande, K. ; Czaja, I. ; Wieneke, U. ; Schell, J., *Science* 1996, 370.
12. S. A. Casson, *The Plant Cell Online* 2002, 14, 1705.
13. N. N. Narita, S. Moore, G. Horiguchi, M. Kubo, T. Demura, H. Fukuda, J. Goodrich, H. Tsukaya, *Plant J* 2004, 38, 699.
14. S. Djakovic, J. Dyachok, M. Burke, M. J. Frank, L. G. Smith, *Development* 2006, 133, 1091.
15. J. Ma, B. Yan, Y. Qu, F. Qin, Y. Yang, X. Hao, J. Yu, Q. Zhao, D. Zhu, G. Ao, *J Cell Biochem* 2008, 105, 136.
16. X. Dong, D. Wang, P. Liu, C. Li, Q. Zhao, D. Zhu, J. Yu, *J Exp Bot* 2013, 64, 2359.
17. T. P. Kondo, S. ; Zanet, J. ; Benrabah, E. ; Valenti, P. ; Hashimoto, Y. ; Kobayashi, S. ; Payre, F. ; Kageyama, Y., *Science* 2010, 329, 336; M. I. Galindo, J. I. Pueyo, S. Fouix, S. A. Bishop, J. P. Couso, *PLoS Biol* 2007, 5, e106; M. M. Kessler, Q. Zeng, S. Hogan, R. Cook, A. J. Morales, G. Cottarel, *Genome Res* 2003, 13, 264; K. Hanada, X. Zhang, J. O. Borevitz, W. H. Li, S. H. Shiu, *Genome Res* 2007, 17, 632.
18. V. Olexiuk, G. Menschaert, *Adv Exp Med Biol* 2016, 926, 49.
19. S. A. M. Slavoff, A. J. ; Schwaid, A. G; Cabili, M. N. ; Ma, J. ; Levin, J. Z. ; Karger, A. D. ; Budnik, B. A. ; Rinn, J. L. ; Saghatelian, A, *nature chemical biology* 2012, 59.
20. J. Ma, C. C. Ward, I. Jungreis, S. A. Slavoff, A. G. Schwaid, J. Neveu, B. A. Budnik, M. Kellis, A. Saghatelian, *J Proteome Res* 2014, 13, 1757.
21. S. Degroeve, L. Martens, *Bioinformatics* 2013, 29, 3199.
22. V. Olexiuk, J. Crappe, S. Verbruggen, K. Verhegen, L. Martens, G. Menschaert, *Nucleic Acids Res* 2016, 44, D324.

MP-12 / PROTEOMIC CHARACTERIZATION OF HUMAN PLATELET LYSATE BY LABEL-FREE NLC-MS/MS: TOWARDS THE EFFECTORS OF ITS BENEFICIAL EFFECT IN REGENERATIVE MEDICINE

Clizia Chinello⁽¹⁾ - *Federica Re*⁽²⁾ - *Camillo Almici*⁽³⁾ - *Andrew James Smith*⁽¹⁾ - *Francesca Paoletta*⁽⁴⁾ - *Gina Lisignoli*⁽⁴⁾ - *Isabella Piga*⁽¹⁾ - *Martina Stella*⁽¹⁾ - *Mariafrancesca Mazza*⁽¹⁾ - *Vanna Denti*⁽¹⁾ - *Luciana Sartore*⁽⁵⁾ - *Domenico Russo*⁽²⁾ - *Fulvio Magni*⁽¹⁾

University of Milano-Bicocca, Department of Medicine and Surgery, Clinical Proteomics and Metabolomics Unit, Vedano al Lambro (1) - Brescia University, ASST Spedali Civili Brescia, Unit of Blood Diseases and Stem Cell Transplantation, DPT of Clinical and Experimental Sciences, Brescia (2) - ASST Spedali Civili of Brescia,

Laboratory for Stem Cells Manipulation and Cryopreservation, Department of Transfusion Medicine, Brescia (3) - IRCCS Istituto Ortopedico Rizzoli, SC Laboratorio di Immunoreumatologia e Rigenerazione Tissutale, Bologna (4) - University of Brescia, Department of Mechanical and Industrial Engineering, Veduggio al Lambro (5)

Keywords: Regenerative medicine, Proteomics, platelet lysate, scaffold, MALDI imaging

Introduction

Hydrogel-forming polymeric scaffolds (HyPS) provides a favourable 3D biomimetic and biodegradable environment for Human Mesenchymal Stromal Cells (hMSCs) [1]. Recently, in HyPS engineered with hMSCs, human platelet lysate (hPL) has been shown as highly effective to stimulate the hMSCs cell-growth and osteo-chondral differentiation [2-4]. The proteins responsible for its beneficial effect have been not yet clarified and will be herein investigated by MS.

Methods

Quali-quantitative evaluation of hPL (n=4) compared to PlateletsPoor Plasma (PPP)(n=4) was achieved by label-free shotgun proteomics. IgG and albumin depleted and not depleted samples were trypsinized [5] and analysed by nLC-UHRTOF. Data was submitted to PEAKS studio and Mascot.

Proteomic images of chondrogenic differentiated hMSCs from bone marrow (BM) and adipose tissue (AT), grown in HyPS for 21 days with hPL, were obtained by RapifleX MALDI Tissue Typer™ [6].

Results

About 450 and 570 protein IDs were identified in not depleted and depleted hPL and PPP, respectively. 59 of them resulted specific of hPL and not present in PPP independently from depletion, representing a possible panel of factors involved in cell stimulating effect of this supplement. Relative quantification highlighted 7 and 43 proteins significantly altered in their abundances comparing hPL vs PPP in not depleted and in depleted samples (fold change>1.5; p<0.05; ≥2 unique peptides). 3 of these differences were in common, including a critical proliferation regulator. The functions and networks of the proteins of interest were investigated and compared. In situ proteome evaluation of AT-hMSCs and BM-hMSCs grown in HyPS in vitro models for 21 days in presence of hPL and chondrogenic medium was also provided. Histology evaluation was also performed with H&E. The molecular images of trypsinized samples were co-registered with the stained counterparts and co-localization of specific peptide signals with differentiated cells were enlightened.

Conclusions

Shotgun analysis and label-free quantification of hPL and PPP allowed to 'proteomically' characterize the supplement and to highlight possible candidates responsible for the beneficial effect onto the proliferation and differentiation of in vitro models. MS imaging analysis of seeded scaffold permitted also to recognize proteomic profiles specific of chondrogenic differentiated cells for in situ monitoring cell expansion and maturation.

Novel Aspect

Proteomic characterization of hPL and MS imaging analysis give important information to clarify and monitor the hMSCs chondrogenic differentiation in in vitro HyPS3D model.

References

1. Dey K., Agnelli S., et al. *Int. J. Polym. Mater. Polym. Biomat.* (2018).
2. Orlandi C, et al. *Exp. Dermatol.* (2018) Mar 31.
3. F. Re et al., Abstract, International Translational and Regenerative Medicine Conference, 25-27 April (2018), Rome
4. F. Re et al., Abstract number: B338, 44th Annual meeting of the European Society for Blood and Marrow Transplantation, 18-21 March (2018), Lisbon.

5. Raimondo F, et al. Mol. Biosyst. (2015) Jun;11(6):1708-16.
6. Galli M, et al. Biochim. Biophys. Acta. Jul;1865(7):817-827 (2017)

Fundings

The research leading to these results has received funding from FAR 2014–2017, Fondazione Gigi & Pupa Ferrari Onlus, UniBs Health & Wealth Project and Fondazione Comunità Bresciana.

MP-13 / VALIDATION OF A LC-MS ASSAY FOR ASSESSMENT OF ALPHA-DEFENSIN IN SYNOVIAL FLUID

Riccardo Dall'Anese⁽¹⁾ - **Paola Pezzati**⁽²⁾ - **Massimo Quercioli**⁽²⁾ - **Fiamma Balboni**⁽³⁾ - **Giovanni Balato**⁽⁴⁾ - **Andrea Baldini**⁽⁵⁾ - **Giuseppe Bartolini**⁽¹⁾

Buzzi Lab, Analytical chemistry, Prato (1) - Azienda Ospedaliero Universitaria Careggi, Centro Regionale di Riferimento SOD Sicurezza e Qualità, Firenze (2) - Istituto Fiorentino Cura ed Assistenza Firenze (iFCA), Laboratorio Analisi Chimico-Cliniche, Firenze (3) - Università Federico II di Napoli, Scuola di Medicina e Chirurgia, Dipartimento di Sanità Pubblica, Napoli (4) - Istituto Fiorentino Cura ed Assistenza Firenze (iFCA), UO Ortopedia, Firenze (5)

Keywords: Synovial fluid, periprosthetic joints infections, alpha defensin, liquid chromatography-mass spectrometry

Introduction

Synovial fluid Alpha defensin is regarded as a potential biomarker for rapid identification of periprosthetic joints infections (PJI). The availability of a method not prone to interferences, with high sensitivity and specificity, could be of great help to clinicians, to make treatment decisions. LC-MS technology is widely used to measure new biomarkers, due to high selectivity and high sensitivity. We developed and validated a LC-MS method to perform synovial fluid alfa-defensin quantification.

Methods

Synovial fluids (SF) samples were sample left over obtained from Patients referring to a large Italian private hospital. A pre-analytical step was performed based on chemical protocol to obtain a peptide from trypsin digestion of protein. The quantification of peptide are carried out by blank synovial fluid (spike solution) and simulant of synovial fluid (simulant [1]) in range from 0.1 to 100 µg/ml. To develop the LC-MS method, two different platforms (LC-QTOF and LC-MS/MS) were used.

Results

The isolation and identification of the marker peptide from alpha defensin was achieved by LC-QTOF analysis by Auto-MS acquisition of data, after tryptic digestion of synovial fluid sample, positive for the presence of alpha defensin (Positive in according to MSIS criteria). The raw data were elaborated with Spectrum Mill Software, to verify the uniqueness of marker peptide of alpha defensin: the peptide LWAFFCC is identify as a specific marker for the presence of alpha defensin. To evaluate the method performance, LOD and LOQ are calculated with 60 replicates of negative samples, in accordance to CLSI EP17-A2:20012 [2]: the best results are 0.09 µg/ml for LOD and 0.29 µg/ml for LOD performed by negative matrix spike and LC-QTOF analysis. For December 2016 to November 2017, 119 synovial fluid samples are tested (from 119 patients, 92 aseptic and 27 positive by MSIS criteria): the results are cut-off value of 1.0 µg/ml, AUC 0.992 (95% C. I.0.98-1.0), sensitivity 100%, specificity 96.3%, PPV 90% and NPV 100%.

Conclusions

Presently, Alfa defensin measure in SF is performed either by laboratory based alpha- defensin immunoassay or by Synovasure “quick test” (lateral flow test kit). Both methods claim high degree of accuracy in the diagnosis of PJI [3]. LC-MS offers the potential to be an alternative to immunoassays and

offers the highest sensitivity and precision for the identification and detection of analytes, although it still requires high skilled professionals and may need a long TAT.

Novel Aspect

The results of the LC-MS assay validation show a robust method for qualitative and quantitative research of alpha defensin in synovial fluid.

References

1. Bortel E. L., Charbonnier B., Heuberger R., *Lubricants*, 3, 664-686 (2015).
2. CLSI EP15-A3:2014 User verification of precision and estimation of bias; approved guidelines. Third Edition 2014.
3. Bossuyt PM., *Scand J Clin Lab Invest Suppl*, 52, 242-246 (2010).

MP-14 / MASS SPECTROMETRY EVALUATION OF CIRCULATING S100A8 AND S100A9 PROTEINS AS PROGNOSIS BIOMARKERS IN SEPSIS

Christelle Dubois⁽¹⁾ - ***Didier Payen***⁽²⁾ - ***Valérie Faivre***⁽²⁾ - ***Stéphanie Simon***⁽³⁾ - ***François Fenaille***⁽⁴⁾ - ***Christophe Junot***⁽¹⁾ - ***Nathalie Morel***⁽³⁾ - ***François Becher***⁽⁴⁾

Service de Pharmacologie et Immunoanalyse (SPI), CEA, INRA, Université Paris Saclay, F-91191 Gif-sur-Yvette cedex (1) - Department of Anesthesiology and Critical Care, Lariboisière Hospital, University of Paris Denis Diderot 7, Paris (2) - Service de Pharmacologie et Immunoanalyse (SPI), Laboratoire d'Etudes et de Recherches en Immunoanalyse (LERI), CEA, INRA, Université Paris Saclay, F-91191 Gif-sur-Yvette cedex (3) - Service de Pharmacologie et Immunoanalyse (SPI), Laboratoire d'Etude du Métabolisme des Médicaments (LEMM), CEA, INRA, Université Paris Saclay, F-91191 Gif-sur-Yvette cedex (4)

Keywords: sepsis, plasma, S100A8/S100A9, biomarker, PRM

Introduction

Sepsis causes one death every 5.2 seconds worldwide [1]. To achieve rapid patient care, prognostic markers are needed. Predictive value of S100A8/9 complex was recently highlighted with an ELISA[2]. However, antibody affinity to particular isoform is questionable. Our study has assessed the prognostic value of plasmatic S100A8 and S100A9 by combining intact protein mass spectrometry for isoform characterization and bottom-up quantification of total isoforms.

Methods

Proteoform characterization includes immunoprecipitation using specific antibodies coated on magnetic beads, acidic elution and intact protein analysis on a Q-Exactive instrument operated in full MS mode. Bottom-up assay: the antibody-free assay consists in a simple pellet and tryptic digestion protocol using 10 µL of plasma. The Q-Exactive was operated in the PRM mode and four specific peptides were monitored with labeled peptides for quantification.

Results

Intact protein analyses following immunoprecipitation were used to characterize the circulating proteoforms of S100A8 and S100A9. In plasma of septic patients, truncated, acetylated and oxidized isoforms were detected.

Peptides common to all detected proteoforms were selected to quantify the total forms. A simple antibody-free sample preparation was developed for targeted analysis of S100A8 and S100A9 signature peptides. Application to clinical cohort analysis requires robust and quantitative assays. Validation in plasma of septic patients demonstrated good reproducibility, linearity and peptide stability.

A cohort of 45 plasma of septic patients was analyzed with the present assay, in parallel to an ELISA [3]. Remarkably, we found a significant difference in S100A8 and S100A9 concentration between survivor and non-survivor patients ($p= 0.009$). Additionally, excellent correlations between MS analysis and ELISA were observed ($r=0.97$), thus demonstrated similar performance for both assays.

Conclusions

Intact mass spectrometry allowed identification of S100A8 and S100A9 isoforms in patients with sepsis and all isoforms were then quantified by a simple antibody-free bottom-up assay. The present method demonstrated high reproducibility and accuracy, in agreement with the latest bioanalytical recommendations [4]. It should help harmonization of S100A8 and S100A9 quantification across laboratories, with the objective of reaching a better sepsis prognosis.

Novel Aspect

Novel validated LC/MS assay for S100A8 and S100A9 quantification in a cohort of patients with sepsis, provide new significant variables for prognosis in sepsis.

References

1. Fleischmann C., Scherag A., Adhikari NK., Hartog CS., Tsaganos T., Schlattmann P., Angus DC., Reinhart K., *Critical Care*, 19, P21 (2015).
2. Gao S., Yang Y., Fu Y., Guo W., Liu G., *American Journal of Emergency Medicine*, 33, 1278-1282(2015).
3. Ikemoto M., Tanaka T., Takai Y., Murayama H., Tanaka K., Fujita M., *Clinical Chemistry*, 49:4, 594-600 (2003).
4. Grant RP., Hoofnagle AN. *Clinical Chemistry*, 60:7, 941-944 (2014).

MP-15 / CHARACTERIZATION OF THE PIF PEPTIDE (PREIMPLANTATION FACTOR) INTERACTOME IN HUMAN SERUM

Anna Fel⁽¹⁾ - ***Paulina Czaplewska***⁽¹⁾ - ***Katarzyna Macur***⁽¹⁾ - ***Marcel Thiel***⁽¹⁾ - ***Krzysztof Łukaszuk***⁽²⁾ - ***Stanisław Ołdziej***⁽¹⁾

University of Gdansk, Intercollegiate Faculty of Biotechnology University of Gdansk and Medical University of Gdansk, Gdansk (1) - Medical University of Gdansk, Department of Obstetrics and Gynecological Nursing, Medical University of Gdansk, Gdansk (2)

Keywords: infertility, preimplantation factor, LC-MS/MS, affinity chromatography, proteomics

Introduction

Nowadays infertility affects a significant number of couples and for many of them ART procedures provide a unique opportunity of achieving pregnancy. The major factor limiting the success of IVF procedures is embryo implantation. Studies showed that PIF peptide can be responsible for tolerances of the immune system of the mother [1]. The aim of the project is to understand the protein interactions of the PIF peptide in human serum to understand molecular mechanisms responsible for implantation.

Methods

For binding proteins from human serum to PIF peptide we used affinity chromatography. PIF peptide was immobilized on solid support and then human serum was applied to microcolumn with immobilized peptide. The eluted proteins were subjected to bottom up analysis by liquid chromatography tandem mass spectrometry. Fragmentation spectra obtained in connection with available bases of amino acid sequences enabled to create a list of proteins that bind to PIF peptide.

Results

The PIF peptide consisting of the 15 aminoacids residues (H-MVRIKPGSANKPSDD-OH) was chemically synthesized with the variants containing cysteinopentaglycine. These variants enable to prepare affinity column for binding proteins. The bottom up analysis with trypsin as a digestion enzyme enable to identify more proteins in a sample. Fragmentation spectra obtained from LC-MS/MS analysis in connection with SwissProt databases enabled to identify several proteins which bind to PIF peptide. These proteins have potential impact on trophoblast invasion, eventually affecting placental development in the peri-implantation period.

Conclusions

Performed methods of affinity chromatography and LC-MS/MS technique are good tools for analysis of PIF interactome in human serum and other fluids. Presented method reveals that PIF peptide interacts with proteins which are responsible for implantation of the embryo and immune tolerance of the mother. However, holistic approach should be performed to understand whole mechanism of the implantation.

Novel Aspect

The information obtained in this project can be used to develop diagnostic and therapeutic methods for monitoring the development of pregnancy and maintaining endangered pregnancies.

References

1. Barnea, E. R., Kirk, D., Ramu, S., Rivnay, B., Roussev, R., & Paidas, M. J. (2012). Journal of Obstetrics & Gynecology, 207, 313-e1 (2012)

MP-16 / MOLECULAR FINGERPRINTING OF THE INTERSTITIUM PROTEOME IN THE VASTUS LATERALIS MUSCLE IN WOMEN WITH FIBROMYALGIA COMPARED TO HEALTHY CONTROLS USING MD-LC-MS/MS

Bijar Ghafouri⁽¹⁾ - *Maria Turkina*⁽²⁾ - *Nikolaos Christidis*⁽³⁾ - *Kaisa Mannerkorp*⁽⁴⁾ - *Eva Kosek*⁽⁵⁾ - *Malin Ernberg*⁽⁶⁾ - *Björn Gerdle*⁽¹⁾

Linköping University, Department of Medical and Health Sciences, Linköping (1) - Linköping University, Department of clinical and experimental medicine, linköping (2) - Karolinska institutet and Scandinavian Center for Orofacial Neurosciences, Department of Dental Medicine, Stockholm (3) - University of Gothenburg, Department of Neuroscience and Physiology, Gothenburg (4) - Karolinska institutet, Department of Clinical Neuroscience, Stockholm (5) - Karolinska institutet, Department of Dental Medicine, Karolinska Institutet, and Scandinavian Center for Orofacial Neurosciences, Stockholm (6)

Keywords: chronic pain, microdialysis, MD-LC-MS/MS, fibromyalgia, intervention

Introduction

Fibromyalgia (FM) (prevalence of 2-4%) is a multifaceted chronic pain condition. The pathophysiological mechanisms behind FM involve both central and peripheral sensitization. Validated objective markers for involved pathophysiological processes in FM are crucial for diagnosis, but are today unfortunately lacking. The aim of this study were to investigate the interstitium muscle proteome in FM and healthy subjects to identify activated molecular pathways before and after progressive resistance exercise.

Methods

30 FM patients and 30 healthy subjects were included. Microdialysis (MD) was performed in the vastus lateralis muscle using commercial catheter (100 kDa). Dialysate samples were collected every 20 minutes during a four hours session at a flow rate of 5 µL/min. The proteome was analyzed by nLC-MS/MS on LTQ Orbitrap Velos Pro (Thermo). Multivariate statistical analysis and bioinformatics were used to compare and describe the proteome profile between the FM and controls.

Results

The proteome fingerprinting analysis of the interstitial muscle identified proteins involved in inflammatory responses, stress response, immune system processes, muscle contraction and nociception. The insertion of MD catheter causes an acute tissue trauma and nociception. The equilibration period for the muscle was two hours. To investigate if the myalgic muscle reacts differently to the insertion of MD catheter compared to healthy muscle, dialysate samples from the trauma period (20-120 minutes) were analyzed. There was a significantly altered protein pattern in FM during the acute tissue trauma compared to the healthy controls. Important proteins whose levels were lower in FM during the tissue trauma were proteins involved in muscle contraction. The expression levels of proteins involved in inflammation, gluconeogenesis and stress responses were higher in FM compared to healthy subjects. Statistical analysis on the effect of the exercise regime on the proteome profile are ongoing.

Conclusions

This study shows that the nociceptive reaction of the muscle due to the MD catheter insertion in FM is different compared to healthy pain free subjects. The findings in this study together with the findings in plasma from the same cohort contribute to increased knowledge about peripheral and systemic molecular mechanisms involved in the nociceptive mechanisms in FM. Taken together our results can contribute to a better understanding of how to improve outcomes of treatments.

Novel Aspect

MD-LC-MS/MS is a promising tool to provide new insights into the interstitial proteome of muscle and to finding altered molecular pathways in fibromyalgia.

MP-17 / DEEP PROTEOMICS OF SAW-SCALED VIPER (ECHIS CARINATUS SOCHUREKI) VENOM BY INTEGRATED TOP-DOWN AND BOTTOM-UP MASS SPECTROMETRIC APPROACHES

Parviz Ghezellou⁽¹⁾ - *Vannuruswamy Garikapati*⁽¹⁾ - *Seyed Mahdi Kazemi*⁽²⁾ - *Kerstin Strupat*⁽³⁾ - *Alireza Ghassempour*⁽²⁾ - *Bernhard Spengler*⁽¹⁾

Justus Liebig Giessen University, Institute of Inorganic and Analytical Chemistry, Giessen (1) - Shahid Beheshti University, Medicinal Plants and Drugs Research Institute, Tehran (2) - Thermo Fisher Scientific (Bremen) GmbH, Thermo Fisher Scientific (Bremen) GmbH, Bremen (3)

Keywords: top-down proteomics, bottom-up proteomics, snake venom

Introduction

Venoms produced by viper snakes contain medically significant toxins that play a key role in the pathophysiology of snakebites disease. The most widely used techniques for venom analysis is LC-MS/MS based bottom-up (BU) proteomics [1]. On the contrary, intact protein measurement with top-down (TD)MS strategies, offers the richest data, covering the almost complete protein sequence and the ability to identify proteoforms and to localize modulations [2-3].

Methods

Snake specimens were captured and milked in their localities around three different geographical regions of Iran. Size-exclusion chromatography (Bio-Rad NGCQuest Plus) was used for protein fractionation prior to LC-MS/MS. Mixed and fractionated venom samples were subjected to TD and BU proteomics analysis using Ultimate 3000 RP-UHPLC coupled to a Fourier transform orbital trapping mass spectrometry (Q Exactive HF-X, Thermo Scientific GmbH, Bremen, Germany).

Results

The venom toxin composition of three geographically different saw-scaled viper samples were primarily

measured as a venom mix through both, direct TD and BU strategies. The results of these approaches provide an initial overview of the venom toxin family on the protein level. The venom mix measurement illustrated that the venom proteomes were dominated by PLA2 and SVMP. Considering the proteome complexity of snake venoms, pre-fractionation is an essential part of venomomics prior to MS analysis. We therefore de-complexed the venom of three populations by size-exclusion chromatography into 6 fractions for each sample. The fractions were then applied to LC-MS/MS in reduced and non-reduced form (TD) and also digested in solution (BU) for proteomics measurements. The data indicated sequence similarities to more than 300 proteins, belonging to 20 protein and peptide families.

Conclusions

Accurate understanding of venom proteome contents can greatly help to improve the treatment of this neglected tropical disease and also drug discovery. Venom research is constantly developing, based on advanced analytical technology. Implementation of multi-dimensional separation techniques and TD proteomics using high-performance mass spectrometry accelerates the identification and characterization of venom components.

Novel Aspect

First comprehensive qualitative and quantitative proteomics study of saw-scaled viper (*Echis carinatus sochureki*) venom from three different geographical regions of Iran

References

1. Calvete J. J., Expert Rev. Proteomics 11, 315–29 (2014).
2. Kelleher N. L., Anal. Chem. 76, 196 A-203 A (2004).
3. Catherman A. D., Skinner O. S., Kelleher N. L., Biochem Biophys Res Commun 4, 683–693 (2014).

MP-18 / DIFFERENTIAL PROTEOME ANALYSIS OF HUMAN NEUROBLASTOMA XENOGRAFT PRIMARY TUMORS AND MATCHED SPONTANEOUS DISTANT METASTASES

Lorena Haene⁽¹⁾ - *Hartmut Schlüter*⁽¹⁾ - *Laura Heikaus*⁽¹⁾ - *Tobias Lange*⁽²⁾

University Medical Center Hamburg-Eppendorf, Institute of Clinical Chemistry and Laboratory Medicine, Hamburg (1) - University Medical Center Hamburg-Eppendorf, Institute of Anatomy and Experimental Medicine, Hamburg (2)

Keywords: spontaneous metastasis, bioluminescence imaging, infrared laser ablation, proteome analysis

Introduction

Metastasis formation is the major cause for cancer-related deaths and the underlying mechanisms remain poorly understood. In this study we describe spontaneous metastasis xenograft mouse models of human neuroblastoma used for unbiased identification of metastasis-related proteins by applying an infrared laser (IR) for sampling primary tumor (PT) and metastatic tissues, followed by mass spectrometric proteome analysis.

Methods

IR aerosol samples were obtained from ovarian and liver metastases, which were indicated by bioluminescence imaging (BLI), and matched subcutaneous primary tumors. Tissue ablation products were tryptically digested using a filter aided sample preparation protocol [1].

Results

Corresponding histology proved the human origin of metastatic lesions. Ovarian metastases were commonly larger than liver metastases indicating differential outgrowth capacities.

Among ~1,700 proteins identified at each of the three sites, 89 proteins were differentially regulated in ovarian metastases while 290 proteins were regulated in liver metastases. There was an overlap of 26 and

10 proteins up- and downregulated at both metastatic sites, respectively, most of which were so far not related to metastasis such as LYPLA2, ACTL8, EIF4B, LGALS7, GFAP, and ELAVL4. Moreover, we established in vitro sublines from PT and metastases and demonstrate differences in cellular protrusions, migratory/invasive potential and glycosylation.

Conclusions

Our study demonstrates that BLI-guided collection in the LAN-1 neuroblastoma model is suitable to identify novel candidate proteins that contribute to spontaneous metastasis formation in vivo. The number of regulated proteins (showing differential expression levels in PT vs. metastatic sites) decreases in advanced metastases. The derived sublines maintain pro-metastatic features in vitro making them highly attractive tools for future functional studies.

Novel Aspect

Summarized, this study identified potential novel putative drivers of metastasis formation and provides novel sublines of metastatic neuroblastoma.

References

1. Wisniewski, J. R. *Methods Enzymol.* 585, 15-27 (2017)

MP-19 / BIOMARKER CANDIDATES FOR ACUTE MYOCARDIAL INFARCTION IN HUMAN SERUM

Min-Jung Kang - Yunseo Seong

Korea Institute of Science and Technology, Molecular Recognition Research Center, Seoul

Keywords: Low abundance proteins; AMI biomarker candidates; Quantitative analysis; UHPLC-LTQ-Orbitrap; ELISA

Introduction

Acute myocardial infarction (AMI) is a chronic disease and major cause of mortality worldwide, together with cancer. AMI is currently diagnosed based on an electrocardiogram and elevated levels of protein biomarkers such as troponin I, troponin T, myoglobin, creatine kinase, and CK-MB [1]. Once a heart attack occurs, surgery should be done within a short period of time, and the prognosis is very poor. If a novel biomarker in early stage is discovered, then it will help to reduce the burden of life for heart disease patients. We used a mixed non-targeted and targeted approach to discover novel biomarker candidates for AMI from the comparison of healthy and patient serum using quantitative analysis.

Methods

Serum samples of healthy controls (n = 22) and AMI patients (n = 14) were pretreated using an albumin and immunoglobulin G depletion kit before trypsin digestion, and desalting. UHPLC-LTQ-Orbitrap analysis was achieved by data dependent analysis followed by a Mascot search. Comparative and shot gun proteomics developed were modified for discover biomarker candidates [2, 3]. Multi reaction monitoring by UPLC-LTQ-Orbitrap and ELISA were used for the quantitative analysis of biomarker candidates from serum samples for validation.

Results

Thirty four peptide biomarker candidates were screened. Multi reaction monitoring showed that alpha-1-acid glycoprotein 1 level was significantly decreased ($P < 0.05$) and those of antithrombin III, clusterin, obscurin, and titin were significantly increased ($P < 0.001$) in AMI patient samples compared to levels in healthy controls. Four significantly decreases peptides, which were unidentified in Mascot Search, were also discovered ($P < 0.05$) as AMI diagnostic biomarker candidates. We quantified the biomarker candidates

further using commercial ELISA kits. ANOVA showed that obscurin and titin concentrations were significantly elevated in AMI patients compared to levels in healthy controls, together with that of troponin I. Obscurin and titin showed the highest decision performance, with sensitivities of 86% and 85%, specificities of 100% and 94%, and diagnostic criteria of 11.4 ng/ml and 36.1 ng/ml, respectively.

Conclusions

A modified comparative proteomics approach provided excellent results in the search for novel biomarker candidates of AMI. A combination of obscurin and titin could be used for the diagnosis of AMI with high diagnostic performance.

Novel Aspect

Combination of new technology, abundant protein depletion; UHPLC-LTQ-Orbitrap; data dependent analysis; Mascot search, was applied to discover AMI biomarker candidates by comparative proteomics. Obscurin and titin are novel biomarker candidates for AMI patient in early stage.

References

1. Babuin, L. and Jaffe, A. S. CMAJ : Canadian Medical Association Journal 173, 1191-1202(2005).
2. Gorg, A., Weiss, W. and Dunn, M. J. Proteomics 4, 3665-3685(2004).
3. Meissner, F. and Mann, M. Nature Immunology 15, 112(2014).

MP-20 / RHEUMATOID FACTOR-CORRELATED PROTEINS IS COMPLEMENTARY BIOMARKER OF RHEUMATOID ARTHRITIS

Hee-Gyoo Kang - Sora Mun - Arum Park - Ae Eun Seok - You-Rim Lee - Hyo-Jin Kim - Yoo-Jin Lee

Department of Senior Healthcare, BK21 Plus Program, Graduate School, Eulji University, Daejeon

Keywords: Biomarker, Proteomics, Rheumatoid factor, Rheumatoid arthritis

Introduction

Rheumatoid factor (RF) is auto-antibody that is involved in autoimmune such as rheumatoid arthritis. RF is diagnostic biomarker of rheumatoid arthritis. However, single RF factor is not sufficient to reflect all of the state of rheumatoid arthritis. To increase efficacy of diagnosis of rheumatoid arthritis, our aim is to discovery RF-correlated proteins that can be used complementarily with RF.

Methods

For proteomics approach, Serums were classified into RF-low (RF<18 IU/ml), RF-intermediate (18<RF<54 IU/ml), and RF-high (RF>54 IU/ml). Nano-liquid chromatography-tandem mass spectrometry is used to identify and quantify global proteins in RF-low and RF-high subjects. Candidate biomarker proteins were validated using Enzyme-linked immunosorbent assay (ELISA) in RF-low, RF-intermediate and RF-high.

Results

Among identified proteins by nano-liquid chromatography-tandem mass spectrometry, 4 proteins, including cholinesterase, histidine-rich glycoprotein (HRG), lipopolysaccharide-binding protein (LBP) and C-reactive protein, were more than 2-fold differentially expressed with statistically significance ($p < 0.05$). Of these 4 proteins, the ELISA results of HRG and LBP was consistent with the results of nano-liquid chromatography-tandem mass spectrometry. Areas under the ROC curve of HRG and LBP were 0.861 and 0.828, respectively. The correlation between RF and HRG was statistically significant ($p = 0.003$), and that with LBP was also statistically significant ($p = 0.044$).

Conclusions

As a results, HRG and LBP are altered in RF-high serum, which is suggested that HRG and LBP can be

complementarily used with RF. In validation discovery set, the efficacy of HRG and LBP proved in patients with rheumatoid arthritis by ELISA.

Novel Aspect

We propose that HRG and LBP can be useful screening markers for rheumatoid arthritis with RF.

References

1. D. Aletaha., T. Neogi., A. J. Silman., J. Funovits., D. T. Felson., C. O. Bingham., N. S. Birnbaum., G. R. Burmester., V. P. Bykerk., M. D. Cohen., *Arthritis & Rheumatology*, 62, 2569-2581(2010).
2. P. Venables., R. Maini., UpToDate 2012, (2006).
3. K. D. Deane., M. K. Demoruelle., L. B. Kelmenson., K. A. Kuhn., J. M. Norris., V. M. Holers., *Clinical rheumatology*, 31, 3-18(2017).
4. I. B. McInnes., G. Schett., *New England Journal of Medicine*, 365, 2205-2219 (2011).

MP-21 / NOVEL BIOMARKER DISCOVERY IN EXTRACELLULAR VESICLES FROM PLASMA FOR COLORECTAL CANCER

Keiko Kasahara⁽¹⁾ - ***Ryohei Narumi***⁽¹⁾ - ***Yoshihiro Shimizu***⁽²⁾ - ***Keiko Masuda***⁽²⁾ - ***Satoshi Nagayama***⁽³⁾ - ***Takeshi Tomonaga***⁽¹⁾ - ***Yoshiharu Sakai***⁽⁴⁾ - ***Jun Adachi***⁽¹⁾

National Institute of Biomedical Innovation, Health and Nutrition, Laboratory of Proteome Research, Osaka (1) - RIKEN center for Biosystems Dynamics Research, Laboratory for Single Cell Mass Spectrometry, Laboratory for Cell-Free Protein Synthesis, Osaka (2) - The Cancer Institute Hospital of Japanese Foundation of Cancer Research, Department of Colorectal Surgery, Tokyo (3) - Kyoto University, Division of Gastrointestinal Surgery, Department of Surgery, Kyoto (4)

Keywords: extracellular vesicles, plasma, biomarker, colorectal cancer

Introduction

Targeted proteomics using SRM has been proved to be highly sensitive and quantitatively accurate. However, it is expensive and time-consuming to perform large-scale SRM analysis. We introduce a systematic and scalable SRM workflow for discovery and validation of novel biomarkers in plasma extracellular vesicles (EVs) targeting more than 400 proteins effectively.

Methods

We developed a three step strategy for selecting biomarker candidates(BCs). Step 1: Primary BCs were selected from three sources; experimental data, literature data and database. Step 2: Targeted proteomics was performed for selecting secondary BCs from primary BCs. Isolated plasma EVs obtained from patients with colorectal cancer(n=59) and healthy volunteers (n=59) were analyzed in Step 1/2. Step 3: We evaluated the performance of secondary BCs in another sample set.

Results

Step 1; 1,300 proteins were identified from plasma EVs as experimental data. 900 and 5400 proteins were extracted from PubMed literature search and ExoCarta database respectively. From these sources, 450 proteins were selected as primary BCs. Step 2; More than 1,000 stable isotope-labeled peptides were synthesized by our original peptide synthesis method, MS-based Quantification By isotope-labeled Cell-free products (MS-QBiC) [1]. We performed SRM analysis by using synthetic peptides to confirm their retention time and select top 5 highest intensity transitions. Next, spike-in plasma EVs isolated from pooled plasma were analyzed to select measurable peptides. Optimized SRM analysis was performed for selecting secondary BCs in individual plasma EVs. Step 2/3 is now still in the investigation stage.

Conclusions

Although the selection of final BCs is still on going, Step1 / 2 are supported by each unique logical method to choose candidates which is the best for purpose. Our three step strategy would be innovative and practical approach for novel biomarker discovery.

Novel Aspect

We proved the feasibility of targeted proteomics using more than 1,000 stable isotope-labeled synthetic peptides for clinical samples.

References

1. Narumi R. et al., Proc Natl Acad Sci U S A. 113(24):E3461-7 (2016)

MP-22 / QUANTITATIVE PROTEOMIC ANALYSIS OF CERVICOVAGINAL FLUID PROTEOME USING A 2D-NLC-ESI-MS/MS FOR EARLY DETECTION OF PRETERM BIRTH

Young Eun Kim - Kwonseong Kim - Dukjin Kang

Center of Bioanalysis, Korea Research Institute of Standards and Science, Daejeon

Keywords: cervicovaginal fluid; preterm birth; quantitative proteomics; biomarkers

Introduction

Preterm birth (PTB) before 37 weeks of pregnancy is one of major causes of perinatal mortality and neonatal morbidity. Despite of previous studies in developing biomarker(s) of PTB [1-3], unfortunately, the valuable protein maker for the early detection of PTB has not been unveiled. In this study, we performed quantitative profiling of cervicovaginal fluid (CVF) proteome for discovery of PTB-specific marker(s) via online 2D-nLC-ESI-MS/MS with isobaric labeling.

Methods

Comparative proteomic analysis of CVF proteome from PTB and control (its normal counterpart) was carried out using online two dimensional nanoflow liquid chromatography-electrospray ionization-tandem mass spectrometry (2D-nLC-ESI-MS/MS) with Isobaric tags for relative and absolute quantitation (iTRAQ)-based labeling.

Results

From 2D-nLC-ESI-MS/MS experiments, a total of 1346 proteins were identified in both CVF proteome from PTB and control and resulted in identifying 30/72 (up-/down-regulated by 1.5-fold change) proteins that were differently expressed in PTB compared to control. Bases on Gene Ontology (GO) enrichment analysis, it revealed that these significantly changed proteins were functionally involved in structural constituent of epidermis and cytoskeleton, serine-type endopeptidase inhibitor activity and oxygen carrier activity. We further evaluated the PTB biomarker candidates in five PTB individuals and seven control individuals, resulting in 55 proteins showed statistically significant difference between PTB and control. Among these, two proteins (transaldolase and junction plakoglobin) have been previously proposed as candidate PTB biomarkers [3] and other proteins are novel PTB biomarker candidates.

Conclusions

In summary, we analyzed the proteomic difference of human CVF between PTB and control group, which allowed the identification of novel CVF protein biomarkers for diagnosis of PTB. Further investigations will focus on the assessment of these biomarkers (specificity and selectivity) in a large cohort of subjects.

Novel Aspect

This study provided novel CVF protein biomarkers that have the potential to be applied in clinical diagnosis of PTB employing quantitative proteomics.

References

1. Goldenberg RL, Goepfert AR, Ramsey PS., Am J Obstet Gynecol. 192, S36-46 (2005)
2. Pereira L, Reddy AP, Jacob T, Thomas A, Schneider KA, Dasari S, Lapidus JA, Lu X, Rodland M, Roberts CT Jr, Gravett MG, Nagalla SR, J Proteome Res., 6, 1269-1276 (2007)
3. Kacerovsky M, Lenco J, Musilova I, Tambor V, Lamont R, Torloni MR, Menon R; PREBIC Biomarker Working Group 2012-2013., ReprodSci, 21, 283-295 (2014)

MP-23 / ANTIBODY-FREE LC-MS METHODS FOR LOW TO SUB NG/ML QUANTIFICATION OF THE SOLUBLE RECEPTOR OF ADVANCED GLYCATION END-PRODUCTS (sRAGE) IN SERUM

Frank Klont⁽¹⁾ – **Simon D. Pouwels**⁽¹⁾ - **Marrit Hadderingh**⁽¹⁾ - **Marc Joosten**⁽¹⁾ - **Jos Hermans**⁽¹⁾ – **Nico C. van de Merbel**⁽¹⁾ - **Péter Horvatovich**⁽¹⁾ - **Nick H. T. Ten Hacken**⁽²⁾ - **Rainer Bischoff**⁽¹⁾

University of Groningen, Department of Analytical Biochemistry, Groningen (1) - University Medical Center Groningen, Department of Pulmonary Diseases, Groningen (2)

Keywords: Biomarker, Immunoaffinity enrichment, Ion-exchange, LC-MS, Quantification

Introduction

The soluble Receptor of Advanced Glycation End-products (sRAGE) is a decoy receptor moderating pro-inflammatory signals in the lungs. For COPD, sRAGE is considered to be a promising biomarker candidate based on data obtained using antibody-based assays. Here we describe two antibody-free approaches for LC-MS-based sRAGE quantification, one using affimers and one using strong cation exchange (SCX).

Methods

A UPLC-MRM/MS method using affimers, a novel class of affinity ligands, and a method employing SCX solid-phase extraction (SPE) were developed to enable sRAGE quantification in the range of 0.2-10 ng/mL using 50 µL of serum. Both methods were validated according to EMA and FDA guidelines and were compared with a previously developed antibody-based LC-MS method for sRAGE. [1]

Results

Efforts to establish an affimer-based sRAGE method led to an easily automatable workflow, which is capable of analyzing hundreds of samples per week thus having considerable potential for clinical applications. SCX-based sample cleanup at pH 10 enabled sRAGE enrichment with good precision owing to the protein's bipolar charge distribution, despite of having an average pI of 7.8. Compared to the antibody-based assay, our affimer method showed good correlation (R-squared of 0.9) and a bias around-25%, whereas the SCX assay showed moderate correlation (R-squared of 0.5) and reported on average 20% higher results. Plausible explanations for the observed differences include different subsets of sRAGE proteoforms being quantified and varying impacts of sRAGE ligands on the enrichment of this protein.

Conclusions

We developed and validated two antibody-free LC-MS methods for sRAGE quantification. Levels reported by these assays are different from each other as well as from those reported by antibody-based assays. Thereby, these results indicate that sRAGE assays should be critically evaluated and adequately characterized in order for this protein to become useful as a biomarker.

Novel Aspect

We present the unique applications of affimers, a novel class of affinity ligands, and strong cation exchange for low abundant protein enrichment.

References

1. Klont F.,Pouwels S. D., Hermans J., et al., Talanta, 182, 414-421 (2018).

MP-24 / PROTEOME ANALYSIS OF HUMAN NEUROMELANIN GRANULES AND NEURONS IN THE CONTEXT OF DEMENTIA WITH LEWY BODIES

Steffen Kösters - Sarah Plum - Caroline May - Kathrin Marcus

Medizinisches Proteom-Center, Ruhr-University Bochum, Bochum

Keywords: Dementia with Lewy bodies (DLB), Neuromelanin (NM) granules, laser microdissection, LC-MS/MS, proteomic study

Introduction

Neuromelanin is a dark insoluble pigment, which forms so called NM granules with lipids and proteins in the substantia nigra (1). It is widely debated if those granules are neuroprotective or neurodegenerative (2,3). For example, there is a loss of dopaminergic neurons containing NM during the course of DLB (4). A proteome analysis of NM granules and neurons could allow deeper insight into the function of these granules in the pathogenesis of neurological disorders like DLB.

Methods

NM granules were enriched from fresh frozen DLB and control tissue by laser microdissection. A tryptic digestion was performed afterwards. The samples were analyzed by liquid chromatography tandem mass spectrometry (LC-MS/MS). MaxQuant was used for the evaluation of the MS/MS spectra. A data synchronization was conducted with Gene Ontology for further analysis.

Results

Mass spectrometric analysis and subsequent evaluation of the data in this pilot study of four DLB and five control samples led to the identification of around 660 proteins. A significant differential study were performed to give a statement about over- and underrepresented proteins. 116 significantly (t-test <0.05) differential proteins were identified: 72 proteins overrepresented and 44 proteins underrepresented in NM granules from DLB cases. Analysis revealed that lysosomal associated proteins are significantly overrepresented and mitochondrial associated proteins are underrepresented in NM granules of DLB compared to healthy control cases. To affirm these results and increase the statistical significance, the study will be extended by 18 DLB and 20 control cases. Additionally, a simultaneously enrichment and comparison of NM granules and neurons will be performed. This could give deeper insight in the function of NM granules in the pathogenesis of DLB.

Conclusions

The results of this pilot study support the hypothesis of the NM granules formation as described in the literature that NM granules are probably lysosome related organelle species. NM forms in the cytoplasm at physiological temperature and pH as a side product of dopamine synthesis and is incorporated by endosomes or lysosomes (5). Furthermore, 536 proteins were identified in all cases, which can be considered as biochemical essential composition of the NM granules.

Novel Aspect

For the first time a protein biochemical characterization of NM granules in the context of DLB was performed.

References

1. Tribl F. et al. *Molecular & Cellular Proteomics*, 4, 945-957 (2005).
2. Paris I. J. et al. *Neurotoxicity Research*, 16, 271-279 (2009).
3. Shima T. et al. *Free Radical Biology and Medicine*, 23, 110-119 (1997).
4. Iseki E. *Neuropathology*, 24, 72-78 (2004).
5. Plum S. et al. *Scientific Reports*, 6, 37139 (2016).

MP-25 / PROTEOMIC PROFILING OF MELANOMA CELL LINES AND TISSUE USING DATA INDEPENDENT ACQUISITION MASS SPECTROMETRY TO PREDICT TREATMENT RESPONSE AND PATIENT SURVIVAL

Christoph Krisp⁽¹⁾ - ***Robert Parker***⁽¹⁾ - ***Dana Pascovici***⁽¹⁾ - ***Nicholas K. Hayward***⁽²⁾ - ***James S. Wilmott***⁽³⁾ - ***Richard A. Scolyer***⁽³⁾ - ***Mark P. Molloy***⁽¹⁾

Australian Proteome Analysis Facility, Macquarie University, Sydney (1) - Oncogenomics Laboratory, QIMR Berghofer Medical Research Institute, Brisbane (2) - Melanoma Institute Australia, The University of Sydney, Sydney (3)

Keywords: Melanoma, SWATH-MS, Targeted Therapy, Inhibitor Resistance, Survival Prediction

Introduction

Melanoma is the most frequent cause of skin cancer-related deaths and is the most common lethal malignancy in young people (<40 yrs). Inhibition of crucial survival pathways, such as the MAPK pathway, has become a routine in melanoma, however, response rates vary and responders cannot be reliably predicted by genotype alone. We conducted mass spectrometric screening to detect cellular processes that might predict MEK inhibition (MEKi) response.

Methods

Ten subcutaneous lymph node metastatic melanoma derived cell lines with known MAPK status and 32 fresh frozen NRASmut or BRAFmut regional lymph node metastatic melanoma specimens were obtained. Lysed cells and tissues were digested and analyzed by Data-Independent Acquisition (SWATH-MS) using a TripleTOF 5600/6600 mass spectrometer. Cell line viability over ten days in presence of 2 μ M of the MEK inhibitor Selumetinib was determined by a PrestoBlue assay.

Results

Across all cell lines, about 2500 proteins were quantifiable. Principal component analysis segregated melanomas based on in vitro MEKi sensitivity, whereas genotype alone did not. In total, 63 proteins were highly correlative with MEKi response. Survival analysis demonstrated significantly better survival ($p=0.01$) for patients with MEKi sensitive (>7.5 years) versus patients with MEKi resistant cell lines (<1.7 years). About 2000 proteins were quantifiable in the tissue cohort. The MEKi proteome phenotype was also a dominant signal in tissues, however, neither this phenotype nor genotype correlated directly with survival. Dividing patients into active and low or inactive melanogenic tumors and combining it with expression level of two plasma membrane proteins, allowed for accurate segregation between patients with good and poor post-biopsy survival. The three marker panel was further validated on a data set retrieved from a publically accessible data base containing RNA sequencing information of 460 melanoma specimen.

Conclusions

Mass spectrometric profiling of patient derived melanoma cell lines using SWATH-MS led to the identification of predictive markers for MEKi therapy response. In a validation cohort using 32 melanoma tissues, SWATH-MS identified the MEKi response as a major driver of patient segregation. Additionally, a three protein marker panel was identified in the tissue cohort predictive for patient post-biopsy survival independent of genotype and MEKi response phenotype.

Novel Aspect

SWATH-MS profiling demonstrated MEKi phenotypes in cell lines and fresh frozen metastatic melanoma and identified a marker panel predicting patient post-biopsy survival.

MP-26 / UNVEILING THE OUTER MEMBRANE VESICLE PROTEOME OF A. SALMONICIDA: A PHYSIOLOGICAL INSIGHT

Tobias Kroniger⁽¹⁾ - ***Andreas Otto***⁽¹⁾ - ***Bernd Köllner***⁽²⁾ - ***Dörte Becher***⁽¹⁾

Institute for Microbiology, University of Greifswald, Greifswald (1) - Institute for Immunology, Friedrich-Loeffler-Institut, Greifswald-Insel Riems (2)

Keywords: Proteomics; Mass spectrometry; *Aeromonas salmonicida*; Subcellular fractionation; Outer membrane vesicles

Introduction

Aquaculture is one of the most rapidly growing food sectors globally [1]. Still, high economic losses occur due to infections, for example by the gram negative bacterium *Aeromonas salmonicida* as this organism causes a septicemic disease called furunculosis in freshwater fish [2]. In this study we utilized mass spectrometry based proteomics to study the physiology of this pathogen with special interest in the outer membrane vesicles that may interact with the host.

Methods

The *A. salmonicida* strain JF2267 was grown in lysogeny broth liquid medium. The bacteria were harvested in the logarithmic and stationary growth phase. Cytoplasmic, inner membrane, outer membrane and extracellular proteins as well as proteins of outer membrane vesicles were enriched by subcellular fractionation and subsequently analyzed in a mass spectrometer. The resulting data were processed with the MaxQuant [3] and Perseus [4] software.

Results

Overall, app. 50 % of the proteins predicted for this *A. salmonicida* strain were identified, indicating a high proteome coverage. For validation of the applied fractionation protocol the predicted cellular location (according to PSORTb [5]) and the subcellular fraction in which the proteins were identified, were compared. This comparison showed a clear enrichment of the proteins according to their prediction. Further, clear alteration between the logarithmic and stationary phase could be observed on proteome level. The changes in the fraction of the outer membrane vesicles are particularly interesting, as more than twice as many proteins were identified in the vesicles of the stationary phase compared to the logarithmic phase.

Conclusions

In summary, in this study a very good proteome coverage could be achieved by application of a subcellular fractionation for *A. salmonicida*. Moreover, as expected, significant alterations were observed in the comparison of the proteome of the logarithmic and stationary phase. Interestingly, the most striking changes were detected in the vesicle fraction. Therefore, in future experiments, we will have a deeper look into the physiological role of the vesicles.

Novel Aspect

For the first time the protein composition of the outer membrane vesicles of *Aeromonas salmonicida* was determined.

References

1. FAO. 2016, The State of World Fisheries and Aquaculture 2016
2. Toranzo A. E., Magariños B., Romalde J L., Aquaculture, 246, 37-61 (2005)
3. Cox J., Mann M., Nature Biotechnology, 26, 1367-1372 (2008)
4. Tyanova S., Temu T., Sinitcyn P., Carlson A., Hein M. Y., Geiger T., Mann M., Cox J., Nature Methods, 13, 731-740 (2016)
5. Yu N. Y., Wagner J. R., Laird M. R., Melli G., Rey S., Lo R., Dao P., Cenk S. S., Ester M., Foster L. J., Brinkman F. S. L., Bioinformatics, 26, 1608-1615 (2010)

MP-27 / QUANTITATION OF MONOCLONAL ANTIBODY INFlixIMAB IN HUMAN PLASMA BY LC-MS/MS USING FAB-SELECTIVE LIMITED PROTEOLYSIS nSMOL TECHNOLOGY

Neil Loftus⁽¹⁾ - *Alan Barnes*⁽¹⁾ - *Aurore Jaffuel*⁽²⁾

Shimadzu, MS/BU, Manchester (1) - Shimadzu France, MS, Paris (2)

Keywords: mAB, nSMOL, LC-MS/MS, infliximab, plasma

Introduction

Infliximab is a chimeric IgG1 kappa monoclonal antibody, used in the treatment of inflammatory diseases targeting tumor necrosis factor-alpha (TNF). However, studies have identified a substantial proportion of patients (between 30-40%) who fail to respond to anti-TNF therapy. In this study, we describe the use of nSMOL (nano-surface and molecular-orientation limited) proteolysis to selectively target the Fab CDR region of the mAb and LC-MS/MS detection.

Methods

Plasma samples (5µL) were initially captured using FG beads and digested with trypsin to selectively target the Fab CDR region of the mAb, followed by LC-MS/MS detection using reversed phase UHPLC separation (Nexera LC, Shimadzu Corporation using a Shimpack GISS 1.9 µm 50x2.1 mm C18).

Results

A novel method was developed that performed selective proteolysis of the Fab region of mABs using nano-surface and molecular-orientation limited (nSMOL) proteolysis (Iwamoto et al. 2014). The approach decreases sample complexity through mAb protein enrichment from plasma and uses proteotypic peptide sequence specificity for quantitative LC-MS/MS analysis. By restricting digestion, tryptic fragments are limited to the Fab region, allowing a selective quantification of target mAb peptides with higher efficiency and selectivity in combination with MRM measurement. Skyline software (MacCoss Lab, University of Washington) was used to perform in-silico protein digestion and predict candidate peptides and MRM transitions. SINSATHYAESVK, YASEMSGIPSR and DILLTQSPAILSVPGER were selected based on response, linearity and accuracy.

Conclusions

Quantitation of the mAb infliximab in plasma was implemented using standard sample preparation and generic peptide chromatographic conditions. Method development using Shimpack GISS column achieved a cycle time of 8 minutes per sample. nSMOL technology in combination with LC-MS/MS can offer a rapid alternative to mAb quantitation by ELISA. Inclusion of an internal standard in nSMOL& LC-MS/MS increases accuracy and robustness not available in ELISA type assays.

Novel Aspect

Using a Fab CDR restricted trypsin digestion method, infliximab was quantified in plasma by LC-MS/MS resulting in an alternative technique compared to ELISA.

References

1. Iwamoto N., Shimada T., Umino Y., Analyst, 139, 576 (2014).

MP-28 / TOP DOWN MASS SPECTROMETRY WITH MULTIPLE MS/MS STRATEGIES TO IDENTIFY AGE-RELATED PROTEOFORM CHANGES IN TEAR FLUID

Daniel Lopez-Ferrer⁽¹⁾ - *Romain Huguet*⁽²⁾ - *David M. Horn*⁽²⁾ - *Peter P.M. Raus*⁽³⁾ - *Vlad Zabrouskov*⁽²⁾ - *Andreas Huhmer*⁽¹⁾ - *Peter Dem Verhaert*⁽⁴⁾

Thermo Fisher Scientific, Vertical Marketing, San Jose (1) - Thermo Fisher Scientific, Product Management, San Jose (2) - Miro Center, Clinic, Geel (3) - ProteoFormiX, R&D, Beerse (4)

Keywords: biofluids, precision-medicine, top-down proteomics, UVPD, Orbitrap

Introduction

Tear fluid is of key importance to maintain the health of the front of the eye and provide clear vision. In addition, it represents a promising body fluid that can help in the diagnosis and prognosis of various eye diseases. Here we explore the use of top-down proteomics to map the proteoform diversity in the tear fluid. The proposed work flow includes an array of top-down fragmentation techniques on a chromatographic time scale that allows for identification of over 50 proteoforms simultaneously in less than 1 hour. We then applied this workflow to a cohort of sixteen individuals with ages ranging from 18 up to 70 years old to identify molecular changes due to age

Methods

Tear samples were purchased from Lee Biosolutions. Extracted proteins were directly transferred to an injection vial and analyzed by LC-MS using an Orbitrap Fusion Tribrid Instrument modified with UVPD. MS/MS acquisition was performed using ETD, ETHcD, HCD and UVPD fragmentations at a resolution of 120,000 at m/z 200. Data analysis was performed using Thermo Scientific Deconvolution 4.0 and ProSight PD nodes in the Proteome Discoverer 2.2 software. Label free analysis was done using Biopharma Finder 3.0. Statistical analysis to identify differentially expressed or processed proteoforms was done using R.

Results

Tear fluid functions as lubrication of the eye surface, contains a cocktail of antimicrobials to prevent bacterial infections, washes away foreign matter that is deposited in the eye during our daytime activities, and keeps the eye surface clean and smooth. Tears have the potential to be a source of biomarkers for diagnosing eye-related diseases. First, it is a fluid that can easily be obtained in a non-invasive manner. Second, it is localized to the organ of interest to the ophthalmologist (the eye) and not prone to the same dilution problems associated with biomarker fluids like plasma. Third, tear fluid composition changes with the health of the eye (as well as with age).

In this work, we have coupled high sensitivity top-down mass spectrometry to enable the characterization of tear proteoforms in less than 1 h with an elegant combination of bioinformatics algorithms. Our approach consists of a simplified two-step workflow. First, samples are then analyzed over a 30 min LC gradient by MS and MS/MS using various fragmentation mechanisms. And then, MSMS are searched against a protein database. Label free quantitation is performed by aligning the raw files in the chromatographic dimension and extracting the areas under the curve for each of the identified molecular species. Intensities were then normalized and a t-test was performed to identify differential events. PCA analysis shows that samples within the age groups are more similar. And among the significant proteins in the dataset, lacritin is the one that changes most with age. This protein is involved in the promotion of basal tearing and which low levels are associated with dry eye syndrome. Overall, this strategy offers a powerful option for discovery and characterization of potential tear biomarkers that could be used for the screening of diseases, both eye related or other.

Conclusions

From our analyses, we have proven that tear fluid could be of great interest for understanding eye related diseases and stratifying individuals. In addition, this work also proves that top-down proteomics in clinical research presents a great opportunity to inform about the actual changes that are happening, because it enables to correlate exact proteoforms to a given phenotype. And opens the door for the development of proteoform specific assays in the clinical setting.

Novel Aspect

Tear proteoform biomarker discovery in hours using label free proteomics and an Orbitrap Fusion Lumos.

MP-29 / ADVANCES IN CLINICAL PROTEOMICS FOR ANALYSIS OF THYROID FINE NEEDLE ASPIRATION BIOPSIES: EVALUATING PROTEOMIC STABILITY IN PRESERVATIVE SOLUTIONS

Isabella Piga⁽¹⁾ - *Giulia Capitoli*⁽²⁾ - *Silvia Tettamanti*⁽¹⁾ - *Vanna Denti*⁽¹⁾ - *Andrew James Smith*⁽¹⁾ - *Clizia Chinello*⁽¹⁾ - *Martina Stella*⁽¹⁾ - *Davide Leni*⁽³⁾ - *Mattia Garancini*⁽⁴⁾ - *Stefania Galimberti*⁽²⁾ - *Fulvio Magni*⁽¹⁾ - *Fabio Pagni*⁽⁵⁾

University of Milano-Bicocca, Department of Medicine and Surgery - Clinical Proteomics and Metabolomics Unit, Vedano al Lambro (1) - University of Milano-Bicocca, Department of Medicine and Surgery - Centre of Biostatistics for Clinical Epidemiology, Monza (2) - San Gerardo Hospital, Department of Radiology, Monza (3) - San Gerardo Hospital, Department of Surgery, Monza (4) - University of Milano-Bicocca, Department of Medicine and Surgery - Section of Pathology, Monza (5)

Keywords: thyroid lesions; Fine needle aspiration biopsies, proteomics, cytological sample stability, mass spectra similarity index

Introduction

Thyroid nodule lesions is one of the diseases diagnosed using liquid-based cytology [1]. Such samples are collected via fine-needle aspiration biopsy (FNAB) and deposited into preservative solutions in order to maintain their morphological integrity. However, the proteomic stability of these preserved samples is yet to be investigated and must be ascertained in order for them to be reliably employed in proteomic studies aimed at biomarker discovery.

Methods

Thyroid FNABs were collected from 14 patients (San Gerardo Hospital, Monza, Italy) and transferred into CytoLyt solution, centrifuged and re-suspended in PreservCyt solution. Cytospin spots have been positioned onto ITO-conductive slides and MALDI-MSI intact proteins analysis was performed using an ultrafleXtreme MALDI-TOF/TOF (Bruker Daltonik). Spectra pre-processing and data analysis were performed using the open-source R software v. 3.4.3.

Results

Each FNAB was split into several samples in order to investigate the experimental repeatability (intra-day and inter-day) of the proteomics analysis and the cytological samples stability after 7, 14 days and 2 months in PreservCyt and after 7 days in CytoLyt at 4°C. All samples were compared with the one prepared at t0. Mass spectra similarity was evaluated by using two score systems. The score S3, derived from a previous study [2], was the sum of three components (fit, retrofit and spearman's correlation). The second score system (S4) ranges 0–4 and includes a fourth feature that measure the overlap, which takes into account the whole shape of the two spectra. Intra-day and inter-day CV were very low within ranges of 8.64%-12.03% for S3 and of 7.37%-10.43% for S4, respectively. The results suggest no substantial deviations from t0 when the cytological samples were stored in PreservCyt until 14 days and in CytoLyt until 7 days. However longer storage time (2 months) in PreservCyt does not preserve the specimens, and the spectra overlap with t0 was only 50%.

Conclusions

This study represents a step forward towards the implementation of MALDI-MSI, combined with a trustworthy and robust sample preparation methodology, into the cytopathology routine, integrating the morphology with the proteomics data to improve the diagnosis. Moreover, this protocol allow simple sample collection and shipment to be used not only for the proteomic MALDI-MSI analysis of thyroid FNABs but also for other biological liquid based specimens.

Novel Aspect

A new similarity score was introduced to equally take into account the number of signals (fit and retrofit)

and their intensities (spearman's correlation and spectra overlap).

References

1. Rossi ED, Zannoni GF, Monceli S, Stigliano E, et al., *Frontiers in Endocrinology*, 3, 1-4 (2012).
2. Hollemeyer K, Altmeyer W, Heinzle E, Pitra C, *Rapid Communication in Mass Spectrometry*, 22, 2751–2767 (2008).

Funding: This work was funded thanks to AIRC (Associazione Italiana per la Ricerca sul Cancro) MFAAG GRANT 2016- Id. 18445.

MP-30 / QUANTIFICATION OF ENDOGENOUS VEGFA LEVELS FROM HUMAN PLASMA USING NANOLC/HIGH RESOLUTION MASS SPECTROMETRY COMBINED WITH IMMUNOCAPTURE ENRICHMENT

Sina Pleiner - Karen Koehler - Joachim Hoke - Sascha Keller

Boehringer Ingelheim Pharma GmbH & Co. KG, DMPK, Biberach a. d. Riss

Keywords: VEGFA, biomarker, immunoprecipitation, LC-HR/MS

Introduction

VEGFA as biomarker for blood vessel formation and tumor progression. Tumor progression is often accompanied by `aberrant angiogenesis`, a process of blood vessel formation with a special permeable structure. Such structures lead to an adequate supply of the tumor with oxygen and nutrients [1]. Vascular endothelial growth factor A (VEGFA) is an angiogenesis stimulating factor and can therefore serve as a prognostic biomarker for several types of cancer [2, 3].

Methods

VEGFA enrichment and analysis by LC-HR/MS. Human plasma samples were mixed with isotopically labelled VEGFA protein as internal standard. Unlabeled, recombinant VEGFA was used to create a calibration curve. Samples were denatured and digested before VEGFA peptides were enriched by an immunocapture step using the KingFisher™ Flex platform. After further SPE clean up, samples were analyzed in the fullscan mode of a Q Exactive™ HF system connected to a nanoLC system (all Thermo Fisher Scientific).

Results

Quantification of endogenous VEGFA levels from human plasma with high precision. We were able to set up and validate a robust LC-MS method for the quantification of VEGFA as a biomarker from human plasma down to 0.7 ng eq/mL (~ 1.8 fmol absolute per sample) with a CV < 10%. Since the method is based on an immunoprecipitation approach using a polyclonal antibody as capturing agent, antibody lot-to-lot variance was examined as well and found to be negligible.

Conclusions

Approach can be applied to patient samples.

The newly developed LC-MS approach can be used to analyze endogenous VEGFA levels and drug treatment profiles of patient samples. As the method is semi-automated, up to 96 samples can be processed in one run.

Novel Aspect

Development of a LC-MS method for VEGFA.

To our knowledge it was the first time that VEGFA was quantified from a complex plasma sample by an LC-MS assay.

References

1. Folkman J., N Engl J Med, 285, 1182-6 (1971).
2. Leung D. W., Cachianes G., Kuang W. J., et al., Science, 246, 1306-9 (1989).
3. Dvorak, H. F., J ClinOncol, 20, 4368-80 (2002).

MP-31 / CHARACTERIZATION OF THE MODIFICATION IN MITOCHONDRIAL PROTEOME DUE TO EIF6 DEPLETION BY SWATH-MS ANALYSIS AND MULTIVARIATE STATISTICAL TOOLS

Elisa Robotti⁽¹⁾ - *Marcello Manfredi*⁽¹⁾ - *Simona Martinotti*⁽¹⁾ - *Daniela Brina*⁽²⁾ - *Alessandra Scagliola*⁽²⁾ - *Elia Ranzato*⁽¹⁾ - *Emilio Marengo*⁽¹⁾ - *Stefano Biffo*⁽³⁾

University of Piemonte Orientale, Department of Sciences and Technological Innovation, Alessandria (1) - Istituto Nazionale Genetica Molecolare "Romeo ed Enrica Invernizzi", Milano (2) - University of Milano, Department of Biosciences, Milano (3)

Keywords: eIF6; mitochondrial proteome; shotgun proteomics; multivariate statistical analysis; biomarker identification

Introduction

Eukaryotic Initiation Factor 6 (eIF6) binds 60S ribosomal subunits, has with an anti-association property, and rate limiting for tumour onset and progression. eIF6 haplo insufficient cells are normal, but not efficiently transformed in vitro. Some lines of evidence have shown that alterations in mitochondria (the main compartments of energy production) contribute to the development of metabolic syndrome.

Methods

The mitochondrial proteome of AML-12 (non-tumourigenic murine liver hepatocytes) cell line, with eIF6 down-regulated by shRNA, and of liver tissues from wild type and +/- mice for eIF6 has been analysed by uHPLC-QTOF-MS/MS exploiting the SWATH-MS acquisition, which is a high throughput label-free method for protein quantitation combining shotgun proteomics with the quantitative accuracy and reproducibility of selected reaction monitoring (SRM).

Results

The data collected were further processed by multivariate statistical tools to identify the effects of eIF6 depletion. Partial least squares-discriminant analysis (PLS-DA) coupled to a variable selection procedure in cross-validation was applied to this purpose, in order to identify the most relevant panel of markers with the best predictive ability. The analysis allowed the identification of panels of markers with a very good predictive ability. Bioinformatic tools applied to the identified markers allowed the identification of deranged pathways.

Conclusions

We found that depletion of eIF6 by shRNA induces profound and varied impact on mitochondrial proteome, impairing the energy production, steering the metabolism toward the up-regulation of aerobic glycolysis and the inhibition of oxidative phosphorylation.

Novel Aspect

We identified interesting deranged pathways in the mitochondrial proteome after eIF6 depletion.

MP-32 / LABEL-FREE PROTEOMIC ANALYSIS OF URINARY EXOSOMES FROM PATIENTS WITH CC-RENAL CELL CARCINOMA: STAGE-SPECIFIC DIFFERENTIAL PROTEIN PROFILE

Lucia Santorelli⁽¹⁾ - *Francesca Raimondo*⁽¹⁾ - *Clizia Chinello*⁽¹⁾ - *Fulvio Magni*⁽¹⁾ - *Marco Grasso*⁽²⁾ - *Giorgio Bovo*⁽³⁾ - *Marina Pitto*⁽¹⁾

University of Milano-Bicocca, Department of Medicine and Surgery, Vedano al Lambro (1) - San Gerardo Hospital, Urology Unit, Monza (2) - Vimercate Hospital, Pathology Unit, Vimercate (3)

Keywords: ccRCC, Urinary Exosomes, Proteomics, Label-free, Biomarkers

Introduction

RCC represents 2-3% of all cancers and its incidence is increasing [1]. Clear cell RCC (cc-RCC) is the most frequent type. Although ccRCC are diagnosed at early stages, their aggressiveness and clinical outcomes result heterogeneous within each staging group; thus, novel predictors are needed. A non-invasive test using urine will have a significant impact on patient management. Urinary Exosomes (UE) are proposed as a useful biomarkers source for ccRCC.

Methods

The aim of this study is to highlight a characteristic UE protein profile of RCC progression. UE pools obtained from patients with cc-RCC classified as pT1 (n=6) and pT3 (n=6), and from age- and sex-matched healthy subjects (n=6) were isolated by ultracentrifugation. The corresponding proteome of each sample was digested using the FASP-based method and analysed by label-free nLC-ESI MS/MS.

Results

The UE proteome of patients with ccRCC (low- and high-stage) was investigated by mass spectrometric approach following correlation of the data with clinical features (TNM) [2]. The MS analysis allowed the identification and quantification of about 500 proteins. Considering only the protein identification, the data show that the UE protein composition remained rather constant in the three groups (pT1, pT3 and CTRLs). However, by evaluating the protein abundance, 230 had a statistically significant variation. The pathway network elaboration revealed that many typical urinary system proteins changed their levels in UE. Furthermore, the results shown that some proteins involved in signaling pathways change their abundance, increasing disease stage. This finding suggest a possible relationship of these quantitative variations with the process of tumour dedifferentiation and spread. Finally, a panel of proteins were selected in order to validate their differential content between ccRCC and Ctrl UE, through immunoblotting.

Conclusions

UE comparative proteomics could represent a promising starting point for the identification of ccRCC biomarkers for prognostic purposes. The data obtained may be a useful finding for the future investigation of the molecular mechanism of ccRCC progression. More important, with the prospective validation in a larger cohort, information derived from this study may lead to the development of ccRCC candidate biomarkers.

Novel Aspect

In this work, for the first time, the proteomic characterization of UE in patients affected by ccRCC at different stages is approached.

References

1. Ana L Teixeira A. L., Dias F., Gomes M., Fernandes M., Medeiros R., J Kidney Cancer and VHL. 2014; 1: 84–98.
2. Delahunt B., Cheville J. C., Martignoni G, et al. The American journal of surgical pathology. 37(10), 1490–1504 (2013).

Acknowledgements

Granted by the “Regione Lombardia” “Spesa sostenuta a valere sul Bando Linea R&S per Aggregazioni cofinanziata dal POR FESR 2014-2020” CUPE47H16001530009 and in part by Fondazione Gigi & Pupa Ferrari Onlus.

MP-33 / MAXIMIZED THROUGHPUT AND ANALYTICAL DEPTH FOR SHOTGUN PROTEOMICS USING PASEF ON A TIMS EQUIPPED QTOF AND A NOVEL LC SYSTEM

Thomas Kosinsky⁽¹⁾ - *Scarlet Koch*⁽¹⁾ - *Markus Lubeck*⁽¹⁾ - *Pierre-Olivier Schmit*⁽²⁾ - *Nicolai Bache*⁽³⁾ - *Ole Bjeld Horning*⁽³⁾ - *Lasse Falkenby*⁽³⁾ - *Heiner Koch*⁽¹⁾ - *Lucy Woods*⁽⁴⁾

Bruker Daltonik GmbH, Proteomics, Bremen (1) - Bruker France, Proteomics, Wissembourg (2) - Evosep, Application, Odense (3) - Bruker Daltonik GmbH, Bremen, Germany (4)

Keywords: timsTOF Pro, PASEF, Evosep One High throughput, Clinical Research

Introduction

The achievable depth in high throughput proteomics analyses is most often limited by the sequencing speed of mass spectrometers, while the nano-LC turnover rate limits the sample throughput. Here we combine a trapped ion mobility spectrometry (TIMS) QTOF capable of Parallel Accumulation Serial Fragmentation (PASEF) acquisition to a front end allowing efficient use of short gradients, the Evosep One to achieve, high depth at high throughput analysis.

Methods

50ng of a peptide digest of a HeLa cell line were separated on an Evosep One (Evosep) and analyzed on a timsTOF Pro instrument (Bruker Daltonics). PASEF acquisition allowed for very high sequencing speed at high sensitivity. HeLa digests were separated on very short gradients of 11.5 and 5.6 min corresponding to 100 and 200 samples/day, respectively. Data were analyzed using PEAKS studio (Bioinformatics Solution Inc.) and MaxQuant (Max-Planck-Institute of Biochemistry).

Results

The use of short nano-LC gradients reinforces the need of a fast MS instrument to guarantee reasonable analytical depth. TIMS provides a further dimension of separation in addition to the retention time and m/z ratio that increases peak capacity, and it releases the ions from the TIMS cell as concentrated ion packets. By applying the PASEF scan mode on those packets we obtain almost 100% duty cycle, increased sequencing speed (> 100 Hz) and high sensitivity leading to a deep coverage of the proteome over 4 orders of magnitude dynamic range. To test if we can profit from the high speed and sensitivity of PASEF we set up very short gradients on the Evosep system with a total runtime of 14.0 min (11.5 min gradient) and 7.1 min (5.6 min gradient), resulting in 100 samples/day and 200 samples/day, respectively. Evotips were loaded with only 50 ng of HeLa tryptic digest and analyzed. Using the 200 sample per day turnover we found more than 1200 protein families and 8000 peptide matches.

Conclusions

Our results demonstrate that PASEF on the timsTOF Pro coupled to the Evosep One enables fast and robust analysis of large sample numbers with potential applications for high throughput screening and clinical research.

Novel Aspect

High Depth, High Throughput Label-free proteomics for clinical research.

MP-34 / TARGETED QUANTITATION OF PROTEINS FOR DISCRIMINATING OBESE FROM NORMAL-WEIGHT ADOLESCENTS BY LIQUID CHROMATOGRAPHY-MASS SPECTROMETRY

Hyunsuk Shin⁽¹⁾ - *Kisoon Dan*⁽¹⁾ - *Sanghoon Song*⁽²⁾ - *Dohyun Han*⁽¹⁾

Seoul National University Hospital, Biomedical Research Institute, Seoul (1) - Seoul National University

Hospital, Laboratory Medicine, Seoul (2)

Keywords: protein quantification, stable isotope-labeled (SIL) peptide, LC-MS/MS, protein biomarker, multiple reaction monitoring (MRM)

Introduction

Targeted quantitation of proteins based on isotope-labeled synthetic peptides with multiple reaction monitoring (MRM) LC-MS/MS has been reported as one of the important technique. The aim of the present study was to select unique peptides of target proteins and to determine the concentration of the target peptides with LC-MS/MS in human plasma to compare between obese and non-obese subjects.

Methods

The process began with a list of target proteins that identified by Tandem mass tag (TMT) strategy in the discovery phase. The next step was to select tryptic peptides and transitions of peptides after collision energy optimization using isotope-labeled peptides with Skyline software. Then, peptides were to analyze the samples by MRM LC-MS/MS method and to integrate and to evaluate the data.

Results

Using detectability test of the endogenous peptides, we selected 100 representative target peptides from the 39 proteins. Stable isotope labeled internal standard peptides for 39 proteins have been synthesized. The concentration of the internal standards for each peptides was spiked to the sample was obtained by level of the unlabeled peptides in the plasma and randomized samples (n=60) were analyzed in triplicate by MRM LC-MS/MS. Relative quantitation of the peptides was carried out by monitoring the ratio of heavy-to-light peptide. Intra-day and inter-day precision of stable isotope-labeled peptides showed less than 15.62%. P value and AUC of the peptides were calculated with Medcal software and AUC of 12 peptides (9 proteins) that represented below 0.05 is showed below 0.0417.

Conclusions

100 Target peptides for the 39 proteins that distinguishing between obese and non-obese subjects were quantified and assessed by MRM LC-MS/MS in human plasma. This method is sensitive, robust and fit for purpose to verification of peptides. Finally, 60 human samples with obese or non-obese were successfully analyzed using our MRM method.

Novel Aspect

Fully integrated MRM LC-MS/MS quantitation approach of targeted proteins to compare obese and non-obese subjects.

References

1. D. C. Liebler, L. J. Zimmerman., *Biochemistry*, 52, 3797-3806 (2013).
2. D. Fenyő, *Methods in Molecular Biology*, 673, 211-222, (2010).
3. M. Rauh., *Journal of Chromatography B*, 883-884, 59-67 (2012).

MP-35 / DEVELOPMENT OF A SANDWICH ELISA FOR THE THROMBIN LIGHT CHAIN IDENTIFIED BY SERUM PROTEOME ANALYSIS

Kazuyuki Sogawa⁽¹⁾ - ***Miho Miyahara***⁽¹⁾ - ***Sakino Yamanaka***⁽¹⁾ - ***Shigetsugu Takano***⁽²⁾ - ***Hideyuki Yoshitomi***⁽²⁾ - ***Masaru Miyazaki***⁽²⁾ - ***Masayuki Ohtsuka***⁽²⁾ - ***Fumio Nomura***⁽³⁾

Azabu university, Biochemistry, Kanagawa (1) - Chiba University, General Surgery, Chiba (2) - Chiba University Hospital, Clinical Mass Spectrometry and Clinical Genetics, Chiba (3)

Keywords: pancreatic cancer, biomarkers, serum, TMT, early detection

Introduction

Pancreatic ductal adenocarcinoma (PDAC) accounts for 85–90% of all pancreatic tumours. The median survival of all PDAC patients is less than 6 months, and the 5-year-survival rate is 3–5%[1]. One of the most crucial reasons for the poor prognosis is the lack of early diagnostic markers for PDAC. To overcome this and improve the outcomes of patients with PDAC, there is an urgent need to identify highly sensitive and specific markers for early detection.

In the present study, we employed the TMT method to generate comparative protein profiles of sera samples obtained from pre- and postoperative PDAC patients. Furthermore, we compared the serum levels of candidate proteins to evaluate the ability to discriminate between PDAC and healthy controls. Serum CA19-9 cannot discriminate the PDAC from other diseases such as pancreatitis and major gastroenterological cancers including biliary tract cancer (BTC).

Methods

Using tandem mass tag labelling and LC-MS/MS, we performed comparative analyses of pre- and postoperative sera from PDAC patients to identify specific serum biomarkers for PDAC. In validation studies, we evaluated the discriminatory power of candidate proteins.

Results

Among the 302 proteins analysed, 20 were identified as potential biomarkers, with C4b-binding protein α -chain (C4BPA) and polymeric immunoglobulin receptor (PIGR) being selected for further analysis. The sera levels of C4BPA and PIGR were significantly higher in the preoperative PDAC patients than in the postoperative ones ($P < 0.008$, $P < 0.036$, respectively). In addition, serum C4BPA levels, but not PIGR, in patients with PDAC were significantly higher than those in healthy controls as well as in patients with pancreatitis and other malignancies including biliary tract cancers (BTC) ($P < 0.001$). The respective area under the receiver operator characteristics (ROC) curve (AUC) was 0.860 for C4BPA, 0.846 for CA19-9 and 0.930 for the combination of C4BPA and CA19-9 in PDAC vs non-cancer individuals. The respective AUC was 0.912 for C4BPA, 0.737 for CA19-9 in Stages I and II of PDAC, 0.854 for C4BPA and 0.264 for CA19-9 in PDAC vs BTC.

Conclusions

We have demonstrated that C4BPA is a novel serum biomarker for detecting early stage PDAC, as well as for distinguishing PDAC from other gastroenterological cancers. Further analysis in large cohort studies will warrant C4BPA as a promising biomarker of PDAC in clinical use.

Novel Aspect

C4BPA is a novel serum biomarker for detecting early stage PDAC.

References

1. Jemal A., Siegel R., Xu J., Ward E., A Cancer Journal for Clinicians, 60, 277–300 (2010).

For information please contact: sogawa@azabu-u.ac.jp

MP-36 / USE OF THE MALDI BIOTYPER SYSTEM WITH MALDI-TOF MS FOR RAPID IDENTIFICATION OF MICROORGANISMS CAUSING BACTERIAL URINARY TRACT INFECTION FROM URINE SAMPLES

Maho Yano ⁽¹⁾ - **Kazuyuki Sogawa** ⁽¹⁾ - **Hiroto Maeda** ⁽²⁾

Azabu University, Biochemistry, Kanagawa (1) - Maeda Veterinary Hospital, Surgery, Hokkaido (2)

Keywords: direct identification, LUTD, MALDI-TOF MS, rapid BACpro, urine

Introduction

With the increasing number of cats maintained as pets, the opportunities to treat cats with lower urinary tract disease have recently increased in the clinical veterinary field. Several studies have reported urinary tract infection of felines with LUTD accounting for 2–10% of cases [1]. We performed direct identification of bacteria in urine using pretreatment kits for the direct application of positive blood culture bottles to MALDI-TOF MS, aiming to improve the low identification rates of *E. faecalis* and 2-morphology colony types.

Methods

Urine samples collected from 50 cats with bacterial cystitis brought to Maeda Veterinary Hospital between August 10, 2015 and March 31, 2017 were used in the study. Sample preparation of the urine was performed using the MALDI Sepsityper kit and rapid BACpro. To identify the isolates, MALDI-TOF MS was performed on the AutoFlex TOF/TOF mass spectrometer.

Results

MALDI-TOF MS seems to require high bacterial counts to be able to provide reliable scores. Experimental inocula with at least 1.0×10^5 CFU/mL were required to obtain reliable scores for *E. coli* by MALDI Sepsityper Kit and rapid BACpro, and the required inocula were even higher for *E. faecalis* by MALDI Sepsityper Kit (5.0×10^6 CFU/mL) and BACpro (1.0×10^5 CFU/mL).

Of the 50 urine specimens, the growth of colonies was observed in 37 specimens: 29 specimens with single colony morphology, 8 specimens with two colony morphology. Thirteen specimens did not grow in culture, and MALDI-TOF MS did not identify any significant protein profile in any of these cases.

Single colony morphology correctly identified 25 isolates (86.2%) using the MALDI Sepsityper Kit and 25 isolates (86.2%) with the rapid BACpro. In 4 cases in which identification was not possible, the bacterial count was 1.0×10^5 CFU/mL or lower, which was below the detection limit of sensitivity. In the specimens with two-colony morphology, the MALDI Sepsityper Kit identified only one species (*Citrobacter freundii*, *Enterobacter cloacae*, *Klebsiella pneumoniae*). The rapid BACpro was capable of identifying 2 species in 7 urine specimens. It was not able to identify any bacterium in one case, in which *E. faecalis* and *E. cloacae* were present, and the *E. cloacae* count was 1.0×10^5 CFU/mL, which was below the detection sensitivity.

Conclusions

The MALDI-TOF MS method employing the rapid BACpro for pretreatment is a quick and reliable method for the identification of bacteria from infected urines, whose shortened analysis time enables an earlier and more accurate selection of antibiotics for feline LUTD treatment.

Novel Aspect

The identification of bacteria by the MALDI-TOF MS method using the rapid BACpro for pretreatment is possible when the bacterial level is 1.0×10^5 CFU/mL or higher.

References

1. Gunn-Moore D., Journal of Feline Medicine and Surgery, 5, 133-138 (2003).

For information please contact: sogawa@azabu-u.ac.jp

MP-37 / TOPOLOGICAL DISTRIBUTION OF LIPID BIOMARKERS IN ADIPOSE TISSUE OF METABOLIC-HEALTHY AND UNHEALTHY SUBJECTS

Roel Tans⁽¹⁾ - **Rieke Bande**⁽¹⁾ - **Arno van Rooij**⁽¹⁾ - **Billy Molloy**⁽²⁾ - **Ron M. A. Heeren**⁽³⁾ - **Rinke Stienstra**⁽⁴⁾ - **Cees Tack**⁽⁴⁾ - **Hans Wessels**⁽¹⁾ - **Alain van Gool**⁽¹⁾

Radboud university medical center, Translational Metabolic Laboratory, Nijmegen (1) - Waters Corporation, Scientific Operations, Wilmslow (2) - Maastricht University, Maastricht MultiModal Molecular Imaging Institute, Maastricht (3) - Radboud university medical center, Department of Internal Medicine, Nijmegen (4)

Keywords: adipose tissue inflammation, oxylipins, imaging mass spectrometry, targeted MRM

Introduction

Adipose tissue inflammation is a well-observed consequence of obesity and a prologue of the metabolic syndrome. It may result from infiltrating immunoinflammatory cells in adipose tissue, which disturbs metabolic homeostasis and initiates low-grade systemic inflammation. Obesity-induced lipid spillover in adipocytes provokes an aberrant adipokine secretory profile, resulting in the release of increased pro-inflammatory mediators such as cytokines and signaling lipid derivatives. Arachidonic acid is the main source of lipid derivatives (*i.e.* prostaglandins, thromboxanes, leukotrienes, resolvins, protectins and epoxyeicosatrienoic acids), that mediate inflammation. Different enzymatic conversions downstream of arachidonic acid, *i.e.* cyclooxygenase (COX), lipoxygenases (LOX) and cytochrome P450 epoxygenases (P450), yield specific lipid derivatives. The correlation between adipose tissue levels of these lipid derivatives and the metabolic healthy or unhealthy phenotype of the obese subjects is not clear [1]. Moreover, it is unclear to what extent circulating lipids reflect the inflammatory state of the diseases adipose tissue [2]. We hypothesize that we can determine the diagnostic value of these molecular lipids as biomarker for adipose tissue inflammation through the identification of aberrant local lipid expressions in diseased adipose tissue and subsequently connecting these to circulating lipid levels. Therefore, this project aims to analyze circulating pro-inflammatory lipids and oxylipins that reflect adipose tissue inflammation. The topological distribution of the molecular lipid profile in adipose tissue will be analyzed with mass spectrometry imaging (MSI), whereas circulating lipids in matched patient plasma samples will be analyzed with LC-multiple reaction monitoring-(MRM) MS.

Methods

Matched adipose tissue and plasma samples of 10 lean- and 10 obese metabolic healthy subjects were included into this study. For MSI experiments, fresh frozen adipose tissue biopsies were cryosectioned at 10 μm thickness slides. Each cryosectioned biopsy was split into three different slides randomly to minimize batch bias effects. A primary MALDI-MSI analysis on human adipose tissue was performed on a RapifleX MALDI TissueTyper system (Bruker Daltonik GmbH, Bremen, Germany) in both positive and negative ionization mode. Representative tissue sections with aberrant lipid expression profiles from the primary screening were subsequently analyzed with high mass resolution imaging on an Orbitrap Elite mass spectrometer (Thermo Fisher Scientific, Bremen, Germany). Data analysis was performed as described previously [3].

A UPLC-ESI-MRM targeted mass spectrometry assay for oxylipins was developed on a Xevo TQ-S triple quadrupole, in negative ESI mode, coupled to an ACQUITY UPLC system using a Waters ACQUITY HSS PFP, 2.1x100 mm, 1.8 μm as analytical column (Waters Corporation, Milford, Massachusetts, USA). COX metabolites were obtained from Cayman Chemicals (Ann Arbor, Michigan, USA) to assess the method's performance, sensitivity and specificity by plasma spike-in assays with mixed mode solid-phase extraction (SPE) using Oasis MAX $\mu\text{Elution}$ plates (Waters Corporation, Milford, Massachusetts, USA). Heavy labeled COX metabolites were used as reference standards in patient samples. Sample preparation and SPE conditions were performed as described previously [4]. Data analysis was performed with TargetLynx V. 4.1 (Waters Corporation, Milford, Massachusetts, USA).

Results

Lean and obese subjects were well defined based on anthropometric parameters (*e.g.* BMI and total body fat mass), biochemical analysis (*e.g.* plasma glucose, insulin and C-reactive protein) and markers of adipose tissue inflammation (*e.g.* IL-6, IL-8, number of macrophages and crown-like structures). Topological distribution of the COX-, LOX- and P450 metabolites in adipose tissue is being determined. Subsequently,

aberrant lipid profiles in adipose tissue will be analyzed with existing metadata of individual subjects to infer the biological relevance of these metabolites. This holistic approach will indicate those relevant oxylipin profiles that are related to adipose tissue inflammation, which are then selected for targeted oxylipin MRM approach in circulating plasma samples.

MRM methods using specific transitions were developed for analysis of fourteen COX metabolites individually (*i.e.* 6 keto-Prostaglandin F1 α , Thromboxane B2, 11 β -Prostaglandin F2 α , Prostaglandin F2 α , 11 β -13,14-dihydro-15-keto PGF2 α , Prostaglandin E2, Prostaglandin D2, 13,14-dihydro-15-keto Prostaglandin E2, 13,14-dihydro-15-keto Prostaglandin D2, 15-keto Prostaglandin E2, Prostaglandin J2, 15-deoxy- Δ 12,14-PGD2, 15-deoxy- Δ 12,14-Prostaglandin J2 and Arachidonic Acid). MRM methods show good linearity across at least two magnitudes of standard ($R^2 > 0.97$). Purification by SPE was further optimized to extract spiked COX metabolites standards (unlabeled and heavy labeled) in human plasma.

We currently perform assay validation of the MRM assay prior to analysis of patient samples. This pilot study of 10 lean and 10 obese subjects will be evaluated first prior to analysis of all lean ($n=35$) and obese ($n=37$) subjects and type 2 diabetic (T2D) patients ($n=33$). The T2D patients will be included as the metabolic unhealthy subjects.

Conclusions

This approach will allow to identify aberrant lipid expressions in diseased tissue and its connection to circulating lipid levels, which indicates the diagnostic value of these molecular lipid profiles as biomarker for adipose tissue inflammation.

Novel Aspect

The unique availability of matched patients samples (*viz.* adipose tissue biopsies and plasma) enables biomarker research for adipose tissue inflammation. This holistic approach that combines complementary mass spectrometry approaches will strengthen personalized biomarker research. This is the first study that combines MALDI-MSI with LC-MRM-MS to analyze molecular characteristics and location of lipid profiles in the adipose tissue and how this relates to the circulating lipid profiles per individual.

References

1. Iyer, A., et al., Inflammatory lipid mediators in adipocyte function and obesity. *Nat Rev Endocrinol*, 2010. 6(2): p. 71-82.
2. Hotamisligil, G. S., Inflammation and metabolic disorders. *Nature*, 2006. 444(7121): p. 860-7.
3. Scupakova, K., et al., Spatial Systems Lipidomics Reveals Nonalcoholic Fatty Liver Disease Heterogeneity. *Anal Chem*, 2018.
4. Katrin Strassburg, B. J. M., Claude Mallet, André Duesterloh, Igor Bendik, Thomas Hankemeier, James Langridge, Rob J. Vreeken, Giuseppe Astarita, Targeted Lipidomics of Oxylipins (Oxygenated Fatty Acids). Waters application note.

MP-38 / LC-MS/MS-BASED CHARACTERIZATION OF GLIOBLASTOMA CELL RESPONSES TO TYPE I INTERFERONS FOR ESTABLISHING DEFECTS IN ANTIVIRAL MECHANISMS

Irina Tarasova⁽¹⁾ - ***Alesya Tereshkova***⁽²⁾ - ***Elizaveta Solovyeva***⁽¹⁾ - ***Julia Bubis***⁽¹⁾ - ***Vladimir Gorshkov***⁽³⁾ - ***Frank Kjeldsen***⁽³⁾ - ***Irina Ilina***⁽⁴⁾ - ***Sergei Moshkovskii***⁽⁴⁾ - ***Peter Chumakov***⁽²⁾ - ***Mikhail Gorshkov***⁽¹⁾

RAS, Talrose Institute for Energy Problems of Chemical Physics, Moscow (1) - RAS, Engelhardt Institute of Molecular Biology, Moscow (2) - University of Southern Denmark, Department of Biochemistry and Molecular Biology, Odense (3) - Orekhovich Institute of Biomedical Chemistry, Moscow (4)

Keywords: glioblastoma, interferon signaling, oncolytic viruses, LC-MS/MS-based proteomics, label free quantitation

Introduction

Therapy of glioblastoma multiforme (GBM) is challenging due to specific property of glioblastoma stem cells (GSC) to penetrate far from the tumor localization, thus, avoiding therapeutic actions. However, experimental therapy using oncolytic viruses have demonstrated ability to kill GSC and cause stable remissions. In this work, we studied specific defects in GBM antiviral mechanisms, enabling success of virus therapy.

Methods

8 GBM cell cultures treated with interferon α, β were analyzed using Orbitrap Fusion Lumos mass spectrometer coupled with UltiMate 3000 nanoflow LC (Thermo Scientific). Protein identification, post-search validation and label free quantification were performed using X!Tandem [1] and MPscore [2], respectively. Statistical analyses of differential protein regulation followed by GO analysis were implemented using python-based libraries and GOrilla [3].

Results

Tests for GBM cells sensitivity to vesicular stomatitis virus Indiana strain (VSV) have revealed 5 GBM cultures developing resistance to virus infection due to pre-treatment with interferon α or β (IFN α , IFN β). Other GBM cultures have either developed resistance after pre-treatment with IFN α or kept sensitivity to the virus regardless of IFN α and IFN β treatment. Based on these results, 5 GBM cultures were combined into the group with intact type I IFN signaling, while the other 3 cultures were defined as cultures with broken (at least partially) type I IFN signaling. Proteomics-based analysis of untreated and IFN treated GBM cells followed by analyses of differentially regulated protein products of IFN stimulated genes (ISGs) and biological process enrichments have revealed that GBM cells with intact IFN signaling respond to IFN α and IFN β treatment by enhanced production of ISGs proteins with antiviral properties. The number of ISGs proteins is considerably less in the group with broken type I IFN signaling.

Conclusions

Results of proteomic analysis of GBM cell responses to IFN α and IFN β treatment supported the results of VSV sensitivity tests. In our data, a molecular signature of significantly upregulated IFIT1/2/3, OAS1/2/3/L, MX1/2 proteins was a marker of developed resistance to VSV infection. We assume inactivation of STAT1 and STAT2, the key components of IFN signaling; can be responsible for the broken antiviral defense in GBM cells.

Novel Aspect

The interactions between GBM cells, GSC and oncolytic viruses as well as the role of IFN signaling for virus therapy are the studies of primary interest, promising novel and effective therapy against glioblastoma.

References

1. R. Craig, R. C. Beavis, *Bioinformatics*. 20 (2004) 1466–1467. doi:10.1093/bioinformatics/bth092.
2. M. V. Ivanov, L. I. Levitsky, A. A. Lobas, T. Panic, Ü. A. Laskay, G. Mitulovic, R. Schmid, M. L. Pridatchenko, Y. O. Tsybin, M. V. Gorshkov, *J. Proteome Res.* 13 (2014) 1911–1920. doi:10.1021/pr401026y.
3. E. Eden, R. Navon, I. Steinfeld, D. Lipson, Z. Yakhini, *BMC Bioinformatics*. 10 (2009) 48. doi:10.1186/1471-2105-10-48.

For information please contact: iatarasova@yandex.ru

MP-39 / EPI TOPE IDENTIFICATION OF GLYCATED HSA AND THE EXTRACELLULAR REGION OF HUMAN RAGE BY SPR AND AFFINITY-MS SPECTROMETRY

Anna Tramarin⁽¹⁾ - **Marina Naldi**⁽¹⁾ - **Loredana Lupu**⁽²⁾ - **Pascal Wiegand**⁽²⁾ - **Edoardo Fabini**⁽¹⁾ - **Genny Degani**⁽³⁾ - **Laura Popolo**⁽⁴⁾ - **Giancarlo Aldini**⁽⁵⁾ - **Michael Przybylski**⁽²⁾ - **Manuela Bartolini**⁽¹⁾

Alma Mater Studiorum, University of Bologna, Department of Pharmacy and Biotechnology, Bologna, Italy (1) - Steinbeis Centre for Biopolymer Analysis & Biomedical Mass Spectrometry, University of Rostock, Department for Proteome Research, Rüsselsheim, Germany (2) - University of Milan, Department of Biosciences, Department of Pharmaceutical Sciences, Milan, Italy (3) - University of Milan, Department of Biosciences, Milan, Italy (4) - University of Milan, Department of Pharmaceutical Sciences, Milan, Italy (5)

Keywords: VC1-RAGE, glycated human serum albumin, affinity-mass spectrometry, epitope identification, surface plasmon resonance

Introduction

In the pathogenesis of diabetic complications, the interaction between the receptor for advanced glycation end products (RAGE) and its ligands AGEs is evident: however AGE structural properties and the interaction mode is still unknown [1]. Aimed at clarifying AGE structural requirements, the study intended to characterize the binding of VC1 (RAGE ectodomain) and commercial glycated human serum albumin (HSAgly) to gain new insights on AGE-RAGE interaction.

Methods

HSAgly and HSA were covalently immobilized onto carboxymethyl dextran surface and the binding with recombinant VC1 was investigated by means of surface plasmon resonance (SPR) [2]. Then, an affinity column was prepared by immobilizing VC1 on CNBr activated sepharose beads for VC1–HSAgly interaction studies by MALDI-ToF mass spectrometry (MS). Epitope extraction MS[3] on HSAgly tryptic digestion was performed to identify HSAgly peptides involved in the binding.

Results

Two surfaces of HSAgly and HSA were prepared via amine coupling chemistry and a SPR-based assay was used to verify the selectivity of the binding between VC1 and HSAgly in comparison with the unglycated isoform. VC1-HSAgly affinity analyses were carried out and a 1:1 binding model was used to extrapolate the equilibrium dissociation constant of the complex (KD). The binding constant (KD = 4.9 μ M) showed a moderate affinity of VC1 for HSAgly and allowed to make some consideration on the strength of such interaction.

Furthermore, a VC1-based affinity column was prepared and VC1-HSAgly interaction was confirmed by MALDI-ToF MS. HSAgly was tryptic digested and analyzed by ESI-Q-ToF and MALDI-ToF MS in order to characterize the glycation sites on the protein scaffold. Epitopes extraction method allowed the isolation and the characterization of HSAgly peptides involved in the interaction. The analysis of the model structure showed that they belong to discontinuous epitope included in domain I and sub-domain IIA of the protein.

Conclusions

The results of this study enabled the assessment of VC1-HSAgly interaction by affinity-MS. VC1 showed to well discriminate between glycated and unglycated HSA confirming that HSAgly acts as VC1 binder. The determination of the affinity binding constant and the structural identification of epitope involved in such interaction provided a reliable characterization of commercial HSAgly as marker for AGE-RAGE interaction studies.

Novel Aspect

The affinity-MS approach for identification of epitopes and affinity quantification of the complex provided key information on the glycation sites involved in RAGE binding.

References

1. Degani G., Altomare A. A., Colzani M., Martino C., Mazzolari A., Fritz G., Vistoli G., Popolo L., Aldini G., *Redox Biology*, 11, 275–285 (2017).
2. Rich R. L., Day Y. S. N., Morton T. A., Myszkowski D. G., *Analytical Biochemistry*, 296, 197–207 (2001).
3. Iuraşcu O., Belaunzar M., Cozma C., Petrusch U., Renner C., Przybylski M., *Journal of The American*

Society for Mass Spectrometry, 27, 1105–1112 (2016).

MP-40 / THE MOLECULAR PATTERN OF PLASMA PROTEINS AND THEIR CORRELATION TO PAIN INTENSITY AND SENSITIVITY IN WOMEN WITH FIBROMYALGIA – A MULTI-CENTER STUDY

Karin Wåhlén⁽¹⁾ - ***Björn Gerdle***⁽¹⁾ - ***Kaisa Mannerkorpi***⁽²⁾ - ***Eva Kosek***⁽³⁾ - ***Malin Ernberg***⁽⁴⁾ - ***Bijar Ghafouri***⁽¹⁾

Linköping University, Department of Medical and Health Sciences, Linköping (1) - University of Gothenburg, Department of Health and Rehabilitation/Physiotherapy, Gothenburg (2) - Karolinska Institutet, Department of Clinical Neuroscience, Stockholm (3) - Karolinska Institutet, Department of Dental Medicine, Stockholm (4)

Keywords: Proteomics, chronic pain, plasma, multivariate statistics, fibromyalgia

Introduction

The complete knowledge of the molecular mechanisms behind fibromyalgia (FM) is still lacking and few studies have analyzed the peripheral molecular pattern in women with FM. The primary aim of this study was to investigate and compare the plasma proteome profile between women with FM and healthy controls (CON). Secondly, to investigate the correlation between plasma protein patterns and pain intensity (Visual analogue scale, VAS) and pressure pain thresholds (PPT).

Methods

Clinical variables and background data were collected and PPTs were recorded using an electronic algometer. Plasma proteome profile from women with FM (n=30) and CON (n=32) was analyzed by 2-DE combined with MS. The 2-DE results were analyzed with Principal Component Analysis (PCA) and Orthogonal Partial Least Squares – Discriminant Analysis (OPLS-DA) models. Correlation between plasma proteins and pain intensity, pain sensitivity was analyzed in FM using OPLS.

Results

As expected significant differences were found between FM and CON in the clinical variables Fibromyalgia impact questionnaire (FIQ), VAS, PPT, hospital and anxiety depression scale (HADS) and pain catastrophizing scale (PCS). The OPLS-DA model, identified 27 plasma proteins as being important for the group separation between FM and CON. In the OPLS model of pain intensity, 12 proteins were multivariately correlated with pain intensity in the FM group; eight of 12 proteins were associated with higher pain intensity score (VAS>50 mm). In the OPLS model of pain sensitivity in FM, 11 proteins multivariately correlated with PPT; three out of 11 proteins were associated with very low pain sensitivity score (<200 kPa).

Conclusions

In this cohort prominent differences exist in plasma protein patterns between FM and CON. Several plasma proteins correlated with pain intensity and sensitivity and these proteins need to be further identified to confirm their role in FM. Combining clinical and proteomic data with multivariate statistics, potential biomarkers could be identified that can lead to improved clinical assessments, treatments and eventually new pharmacological therapies for FM.

Novel Aspect

Plasma proteomic generate large data set, which in combination with multivariate statistics, is an important tool for identifying potential biomarkers associated with FM.

MP-41 / THERAPEUTIC PROTEINS IN SPORTS DRUG TESTING: DETECTION OF DIFFERENT ACTRII INHIBITORS IN SERUM BY MEANS OF (IMMUNO-)AFFINITY PURIFICATION, TRYPTIC DIGESTION, AND LC-HRMS

Katja Walpurgis⁽¹⁾ - ***Andreas Thomas***⁽¹⁾ - ***Christian Reichel***⁽²⁾ - ***Frank Dellanna***⁽³⁾ - ***Mario Thevis***⁽¹⁾

German Sport University Cologne, Institute of Biochemistry/Center for Preventive Doping Research, Cologne (1) - AIT Seibersdorf Laboratories, Doping Control Laboratory, Seibersdorf (2) - Fresenius University of Applied Sciences, Fresenius University of Applied Sciences, Cologne (3)

Keywords: Doping, ActRII-Fc fusion proteins, Therapeutic antibodies, Affinity purification, LC-HRMS

Introduction

Therapeutic proteins represent an emerging class of pharmaceuticals comprising several drug candidates with potential performance-enhancing properties. Especially inhibitors of the ActRII signaling pathways possess the potential for being misused as doping agent in sports as they were found to have positive effects on muscle mass, bone formation, and red blood cell production[1, 2].

Methods

Ammonium sulfate precipitation, (immuno-)affinity purification, tryptic digestion, and LC-HRMS were employed to develop targeted detection assays for the therapeutic antibody Bimagrumab and the ActRII-Fc fusion proteins Sotatercept (ActRIIA-Fc) and Luspatercept (modified ActRIIB-Fc) in doping control serum samples.

Results

Three complementary assays for the detection of different ActRII inhibitors were successfully developed by using different proteomic techniques [3, 4]. The methods were comprehensively characterized and found to be fit for purpose for doping control routine analysis. In addition to LC-HRMS, ion mobility was employed to facilitate the unambiguous identification of the target peptides. For Bimagrumab, clinical samples were analyzed to demonstrate the applicability of the approach to authentic serum specimens.

Conclusions

Mass spectrometric assays for the detection of emerging protein drugs are of great value for both the pharmaceutical industry and preventive doping research. The presented assays will expand the range of available tests and can readily be modified to include other therapeutic fusion proteins and antibody-based drugs.

Novel Aspect

Although no ActRII inhibitors have obtained clinical approval yet, LC-HRMS detection. Methods for three candidates were proactively developed.

References

1. Fedoruk M. N., Rupert J. L., Scandinavian journal of medicine & science in sports,18, 123-131 (2008).
2. Tsuchida K., Sunada Y., Noji S., Murakami T., Uezumi A., Nakatani M., Mini reviews in medicinal chemistry, 6, 1255-1261 (2006).
3. Walpurgis K., Thomas A., Vogel M., Reichel C., Geyer H., Schänzer W., Thevis M., Drug testing and analysis, 8, 1152-1161 (2016).
4. Walpurgis K., Thomas A., Dellana F., Schänzer W., Thevis M., Proteomics clinical applications, doi: 10.1002/prca.201700120 (2017).

MP-42 / A NOVEL HIGH-THROUGHPUT ANTIBODY-FREE LC-MS/MS ASSAY TO QUANTITATE CHROMOGRANIN-A IN SERUM

Darren Weber⁽¹⁾ - Scott Goldman⁽¹⁾ - Michael McPhaul⁽²⁾ - Michael Caulfield⁽²⁾ - Nigel Clarke⁽¹⁾

Quest Diagnostics, R&D Mass Spectrometry, San Juan Capistrano (1) - Quest Diagnostics, Quest Diagnostics Nichols Institute, San Juan Capistrano (2)

Keywords: Tumor biomarker; diagnostics; clinical mass spectrometry; proteomics

Introduction

Chromogranin-A (CgA) is a protein produced in the secretory granules of neuroendocrine tissue, with increased plasma levels being indicative of neuroendocrine tumors (NET)[1]. Currently, CgA is measured using various immunoassays [2]. However, limitations arising from non-specific binding and a reduced dynamic range pose hurdles in implementing such tests in clinical laboratories. Here we describe a novel, antibody-free LC-MS/MS assay to quantitate CgA in serum.

Methods

Exploiting the acidic properties of CgA, 100 μ L of serum is extracted using an anion exchange solid-phase extraction plate followed by addition of internal standard (IS). The extracted sample is concentrated, enzymatically digested, and a unique peptide to CgA is chromatographically resolved and analyzed by SRM on a Sciex 6500+ QTrap mass spectrometer. The ratio of the analyte peak area to the isotopically labeled IS peak area is used to achieve quantitation.

Results

Chromogranin-A showed linearity across a wide range (50-50,000 ng/mL, $R^2=0.99$), as well as inter- and intra-assay reproducibility (CV 5.2-15.7%). Analytical sensitivity and limit of quantitation were determined to be 35.5 and 50 ng/mL, respectively. This method has no cross reactivity with interfering substances, and shows sample stability up to 5 days at room temperature, up to 14 days at 4C, and up to 31 days at -20C. Samples from 308 patients were analyzed to compare CgA serum values measured by the Cisbio CGA-ELISA-US immunoassay with values from our LC-MS/MS assay. After natural logarithm-transformation of the measurements, a normal distribution was seen in the 308 patients. A paired t-test was performed on these data and the concordance between the assays was measured by the Pearson correlation coefficient and Passing & Bablok curve fitting ($R^2=0.76$; $-0.25 + 1.03x$).

Conclusions

We have developed and validated a fully automated LC-MS/MS assay for quantitating chromogranin-A out of serum. Unlike immunoassays, our workflow is not dependent on capture and detection antibodies, and therefore structural changes in CgA due to denaturation, aggregation, or post-translational processing will have no impact on detection.

Novel Aspect

This is a first-of-its-kind LC-MS/MS assay to detect and quantitate chromogranin-A out of serum.

References

1. Donica H, Malecha-Jędraszek A, Strosławska E, Burska A, et al. *Folia Histochem Cytobiol* 2010;48:603-610
2. van der Knaap RHP, Kwekkeboom DJ, Ramakers CRB, et al. *Pract Lab Med*. 2015;1:5-11.

MP-43 / EXOSOMES PURIFICATION FOR NEW BIOMARKERS DISCOVERY IN CANCER

Federica Anastasi⁽¹⁾ - Marco Cecchini⁽¹⁾ - Liam McDonnell⁽²⁾

Laboratorio NEST, Scuola Normale Superiore, Pisa, Italy (1) - Fondazione Pisana per la Scienza ONLUS, FPS,

Pisa, Italy (2)

Keywords: exosomes, proteomics, nLC-MS/MS

Introduction

Exosomes are cell-derived nanovesicles that are secreted from all types of cells and are present in all body fluids. Tumor derived exosomes contain molecules derived from the tumor, and thus represent an opportunity for health monitoring and diagnosis[1].

Here an exosome proteome purification procedure has been developed with the aim to purify exosomes from low sample amount and to study the correlation between exosome cargo and cancer progression.

Methods

Exosomes were purified from 100 µl serum. Serum was previously depleted from bigger vesicles.

Exosomes were then isolated by precipitation and the resulting pellet was ultra-filtrated on 3kDa spin filters, vesicles were lysed, and proteins concentrated. After protein digestion, peptides were analyzed by bottom-up-LC-MS/MS. To reduce the serum proteins level size exclusion chromatography (SEC) has been investigated as an additional purification step.

Results

350 protein groups were identified, and which included 14 of the top 25 frequently identified exosomal proteins defined by the ExoCarta database. Functional enrichment analysis of the proteomics data was performed for cellular component geneontology. Gene ontology analysis performed by STRING reported the presence of exosomal proteins (194 proteins in GO:0070062, FDR: 1.87e-122), PANTHER reported vesicles as the main category present in the sample (135 proteins). Serum albumin, immunoglobulins and protein of the complement system were also identified in the EVs-UF-precipitated samples, indicating that the samples also contained soluble proteins as impurity. Serum exosomes purified by SEC reported lower level of albumin, immunoglobulins and complement system proteins showing the ability of the columns to further purify exosomes from contaminants.

Conclusions

Our exosomes purification procedure will be applied to purify exosomes from cancer patient serum with the aim to study the role of these vesicles in inter-cellular communication and to study the correlation of exosome cargo and tumor progression. Further scaling down of the procedure will be evaluated so as to be able to analyze an even lower amount of starting sample.

Novel Aspect

This procedure enable the in depth proteome analysis of vesicle cargo starting from very low amount of sample. Proteomics analysis of exosomes y bottom-up-LC-MS/MS have been overwhelmingly performed on exosomes purified by ultracentrifugation, but which requires a lot of sample volume available and thus is not feasible on the serum from individual mice (circa 100 µl available from a single mouse) for studying cancer mouse models [2]. The exosomes purification method reported here enables the analysis of exosomes from individual mice.

References

1. Young H. S., Shane F., Zhang V., and Jun C., *Cancers*, 9:8, Jan 2017
2. Mingrui An et al. *J. Proteome Res.*, 2017, 16 (4), pp 1763–1772

MP-44 / DEVELOPMENT OF A NOVEL LC CONCEPT FOR CLINICAL PROTEOMICS

Nicolai Bache⁽¹⁾ - **Philipp E. Geyer**⁽²⁾ - **Christian Ravnsborg**⁽¹⁾ - **Ole Hoerning**⁽¹⁾ - **Lasse Falkenby**⁽¹⁾ - **Peter Treit**⁽²⁾ - **Sophia Doll**⁽²⁾ - **Igor Paron**⁽²⁾ - **Florian Meier**⁽²⁾ - ⁽¹⁾ - **Matthias Mann**⁽²⁾

Evosep, Odense, Denmark (1) - Max-Planck Institute of Biochemistry, Proteomics and Signal Transduction, Martinsried, Germany (2)

Mass spectrometry based proteomics and metabolomics are fast growing and powerful technologies, with the potential to revolutionize health care and precision medicine. However, separation technologies has so far limited throughput and robustness and thereby prevented omics technologies from being fully integrated and routinely used in a clinical setting. Here we describe a conceptually novel low-flow chromatography system that delivers the robustness and throughput required for clinical applications while maintaining the sensitivity of current nano-flow LC instrumentation. The new system uses four low-pressure pumps in parallel to elute samples from a disposable trap column while simultaneously creating a final but offset chromatographic gradient with embedded analytes.

Initially, two low-pressure pumps form a primary gradient that flows through a disposable trap column thereby eluting analytes of interest while two additional low-pressure pumps positioned right after the trap column are modifying the eluent to create an “offset” that increases the retention of the now embedded analytes at the separation column. The offset gradient with the embedded analytes are moved into a long, narrow holding loop that subsequently is switched in-line with a single high-pressure pump and a separation column. Thus elution from the disposable trap column and gradient formation become decoupled from the high-pressure separation.

We have characterized the performance of the new system regarding cross contaminations (<0.05%, total TIC), retention time shifts and peak widths (<8.5 sec) in over 1,500 HeLa runs. The extremely short overhead time of approximately 3 min/gradient allowed us to measure 60 samples per day (21 min gradient, 3 min overhead time) compared to 40 samples with a standard nano-LC systems (22 min gradient, 15 min overhead time). To benchmark reproducibility, we measured 100 human plasma samples in 40 hours and each resulted in several hundred quantified plasma proteins, including more than 50 FDA-approved biomarkers. Assured of a high reproducibility, we then aimed to find significant biological information in the human plasma proteome. For this purpose, we analyzed a longitudinal weight loss study that included a total of 319 plasma samples, where we found strongest effects on proteins of the lipid metabolism and inflammation system.

MP-45 / ANALYSIS AND CHARACTERIZATION OF GASTROINTESTINAL METABOLISM OF MAJOR SOY ALLERGENS GLYCININ AND β -CONGLYCININ IN AN IN-VITRO MODEL BY LC-MS/MS

Julia Bräcker⁽¹⁾ - *Natalie Lok*⁽²⁾ - *Jens Brockmeyer*⁽¹⁾

University of Stuttgart, Food Chemistry, Stuttgart, Germany (1) - University of Hohenheim, Institute of Food Chemistry, Stuttgart, Germany (2)

Keywords: food allergens, in-vitro digestion model, LC-MS/MS, label free quantification, proteomics

Introduction

Soy is one of the eight major allergenic foods, that can cause severe reactions. But little is known so far about the fate of allergenic soy proteins during gastrointestinal digestion. Beside already published potential epitopic regions of different soy allergens with the ability to bind Immunoglobulin E and cause an allergic reaction, there is still no profound knowledge about kinetics of allergen degradation and characterization of degradation products [1].

Methods

The three phases of human gastrointestinal digestion of dried soy beans were simulated with a standardized in-vitro digestion model [2]. While time dependent degradation of intact soy proteins was analyzed by SDS-

PAGE, the identification of the highly complex set of degradation products was performed by micro-LC-MS/MS on a Q-TOF instrument in a software assisted proteomics approach. Degradation rates were estimated by label-free quantification for each phase.

Results

For major soy allergens β -Conglycinin (Gly m 5) and Glycinin (Gly m 6) degradation was observed during gastrointestinal digestion. Gly m 5 β -subunit showed the highest stability during gastric digestion whereas α -subunits of Gly m 5 and Gly m 6 subunits were largely degraded during 30 minutes of gastric digestion. For both phases highly complex mixtures of low-molecular products (Gly m 5: 121 peptides, Gly m 6: 405 peptides) could be identified. Gastrointestinal digestion seems to follow a multi-step mechanism in repetitive patterns regarding time dependent formation and degradation kinetics of peptides with different length. Investigation of digestion kinetics by label-free quantification verifies that some proteins are degraded into larger peptides at an early digestion state, whereas others are degraded as late as in the end of gastric phase. Regions in the amino acid sequence overlapping with previously identified epitopes could be determined, suggesting that these peptides can trigger allergic reactions after intestinal absorption.

Conclusions

By simulation of gastrointestinal digestion of soy in an in-vitro model time dependent degradation of major soy allergens β -Conglycinin und Glycinin could be determined as well as differences in formation kinetics of degradation products. Additionally, degradation products of both allergens with intact epitopic regions could be identified, keeping the ability of IgE-binding to potentially cause an allergic reaction after absorption by the intestinal mucosa.

Novel Aspect

Time dependent degradation of allergens and formation of potential allergenic products could be observed by simulated gastrointestinal digestion in an in-vitro model.

References

1. Saeed H., Gagnon C., Cober E., Gleddie S., Molecular immunology, 70, 125 (2016).
2. Minekus M, Alminger M., Alvito P., Ballance S., Bohn T., Bourlieu C., Carrière F., Boutrou R., Corredig M., Dupont D. et al., Food & function, 5, 1113 (2014).

MP-46 / PROTEOMIC PROFILING OF EXTRACELLULAR VESICLES DERIVED FROM MESENCHYMAL STROMAL CELLS

Jessica Brandi⁽¹⁾ - **Annalisa Adamo**⁽²⁾ - **Marcello Manfredi**^(3,4) - **Elisa Robotti**^(3,4) - **Emilio Marengo**^(3,4) – **Daniela Cecconi**⁽¹⁾ - **Mauro Krampere**⁽²⁾

Department of Biotechnology, University of Verona, Verona, Italy (1) - Department of Medicine, University of Verona, Verona, Italy (2) - ISALIT, Spin-off of Department of Sciences and Technological Innovation, University of Piemonte Orientale, Alessandria, Italy (3) - Department of Sciences and Technological Innovation, University of Piemonte Orientale, Alessandria, Italy (4)

Keywords: Mesenchymal stromal cells, Extracellular vesicles, Shotgun Proteomics

Introduction

Mesenchymal stromal cells (MSCs) are adult multipotent stem cells used for regenerative medicine and immunomodulation [1]. Extracellular vesicles (EVs) released by MSCs act as their paracrine effectors by delivering transcription factors, oncogenes, mRNA, microRNA, and proteins to target cells [2]. In this study, we have investigated the immunomodulatory capabilities of EVs derived from primed-MSCs (pEVs) compared to control (cEVs), through a SWATH-MS approach.

Methods

Shotgun mass spectrometry (MS) analysis has been performed on pEVs and cEVs in order to characterize and compare their proteomic profile. Proteins differentially expressed were evaluated in silico using different bioinformatic prediction tools. Multivariate pattern recognition methods like Principal Component Analysis (PCA) and Partial Least Squares Discriminant Analysis (PLS-DA) were also employed. Selected proteins and pathways were validated using Western Blot and flow cytometry analysis respectively.

Results

Shotgun proteomic analysis revealed in pEVs, a strong modulation of proteins involved in many processes related to immune response. The results obtained through SWATH-MS analysis, were validated through western blotting, by selection of four candidate proteins according to their immunomodulatory potential. Proteins differentially expressed in p-EVs compared to c-EVs were mapped to terms in the GO database. Biological process classification showed enriched terms categories such as “leukocyte migration”, and “leukocyte cell-cell adhesion” suggesting a possible involvement of pEVs in modulating immune effector cells activity. Enrichment analysis of pEVs downregulated proteins revealed that the top four identified pathways are: “PI3K-AKT signaling pathway”, “regulation of actin cytoskeleton”, “focal adhesion” and “leukocyte trans-endothelial migration”. The expression and phosphorylation profile of different components of PI3K-AKT signaling pathway was evaluated by flow cytometry in activated B lymphocytes. Our results showed that, the treatment with pEVs, induced a significant down regulation of the pathway compared to cEVs.

Conclusions

The molecular characterization at protein level of pEVs derived from MSCs, and the consequent functional validation of candidate pathways in activated B lymphocytes, will suggest novel therapeutic perspectives for the treatment of inflammatory and autoimmune diseases.

References

1. Cell Death and Disease (2016) 7, e2062
2. Cell Stem Cell. 2009 July 2; 5(1): 54–63

MP-47 / DIFFERENTIATION OF HYPERTENSIVE DISORDERS DURING PREGNANCY BASED ON URINE PEPTIDOME PROFILING BY HIGH RESOLUTION MASS SPECTROMETRY

Anna Bugrova⁽¹⁾ - ***Alexey Kononikhin***⁽²⁾ - ***Victoria Sergeeva***⁽¹⁾ - ***Maria Indeykina***⁽¹⁾ - ***Igor Popov***⁽²⁾ - ***Natalia Starodubtseva***⁽³⁾ - ***Vladimir Frankevich***⁽³⁾ - ***Evgeny Nikolaev***⁽⁴⁾

Emanuel Institute for Biochemical Physics, Russian Academy of Sciences, Moscow, Russian Federation (1) - Moscow Institute of Physics and Technology, Dolgoprudny, Moscow, Russian Federation (2) - V. I. Kulakov Research Center for Obstetrics, Gynecology and Perinatology, Ministry of Healthcare of the Russian Federation, Moscow, Russian Federation (3) - Skolkovo Institute of Science and Technology, Skolkovo, Moscow, Russian Federation (4)

Keywords: peptidomics, urine, preeclampsia

Introduction

Accurate differential diagnosis of hypertensive disorders is important during pregnancy, because it affects the outcome for the mother and fetus. The high potential of urine peptidomics has been demonstrated for clinical purposes. Several MS based approaches have been successfully applied for urine peptidome analysis and potential biomarker search. Herein we apply LC MS to identify a peptide panel for differentiation of hypertensive disorders during pregnancy.

Methods

A case-control study was designed for 64 women from 4 groups: preeclampsia (PE), superimposed PE, chronic hypertension, and healthy pregnant individuals. The urine peptidome fraction was separated using ultracentrifugation, purified by gel-filtration, and analyzed with nano-HPLC coupled to a high resolution 7T LTQ-FT Ultra (Thermo) mass-spectrometer and semiquantitative label free analysis was performed.

Results

Analysis of the LC-MSMS data from different groups of patients showed that 36 peptides were common for them all. 22 belonged to alpha-1-chain of collagen I, 9 were from alpha-1-chain of collagen III, 2 from alpha-2-chain of collagen I, 1 from alpha-1/2-chain of collagen I, 1 from alpha-1-chain of collagen I/XVIII and one from uromodulin. Patients with hypertensive disorders had 34 common peptides: 12 from alpha-1-chain of collagen I, 10 from fibrinogen alpha-chain, 8 from alpha-1-chain of collagen III, and 4 per other types of collagen. The most important is a group of 16 peptides common for PE and superimposed PE cases. It consists of peptides from alpha-1-chain of collagen I, alpha-2-HS-glycoprotein, apolipoprotein A-I – one per each and the rest 13 belong to alpha-1-antitrypsin. Fibrinogen peptides in the urine point to inflammation processes and those of alpha-1-antitrypsin – to changes in metabolic ways of protein degradation.

Conclusions

Fragments of alpha-1-antitrypsin confirmed their importance as markers of pre-eclampsia suggested earlier. Statistical data analysis by non-parametric Kruskal-Wallis and Mann-Whitney tests showed the presence of 12 peptides differentiating at least one pair of groups.

Novel Aspect

Comparative analysis revealed 12 peptides, which could be used as a diagnostic panel for confident discrimination of various hypertensive disorders.

MP-48 / COMPARATIVE ANALYSIS OF LECTIN BASED GLYCOPROTEINS BETWEEN YOGA GROUPS

J. Kang - Y. Jung - Wonryeon Cho

Wonkwang University, Department of Chemistry, Iksan, Republic of Korea

Keywords: Elderly yoga, Healthy plasma, Comparative glycoproteomics, Lectin affinity chromatography

Introduction

As the elderly population grows, yoga gets more popular as one of the exercises to improve the quality of life and health of the elderly. Glycoproteins in human plasma are known to be associated with multiple diseases, cell differentiation, aging, and etc. This study intends to examine the efficacy of yoga by identifying and comparing the glycoproteins in human plasma from elderly yoga groups before and after participation in a fall-prevention yoga program.

Methods

Human plasma were collected from both Advanced level group and Beginner level group before and after participation in a fall-prevention yoga program. The plasma were pooled and affinity-selected using self-packed LEL columns. The captured glycoproteins from each plasma were then trypsin digested and deglycosylated with PNGase F. After that, proteins were identified with nLC-MS/MS and database searches in both groups.

Results

Plasma glycoproteins were identified from both Beginner level group and Advanced level group before and after 4-month of yoga practice. The average age and sex ratio were equally matched in each group for the

same experimental conditions.

The results show that the number of plasma glycoproteins in Advanced level group is less than that of Beginner's group, demonstrating an inverse correlation between number of plasma glycoproteins and yoga experience.

Conclusions

As the result of 4-Monthyoga Program, the numbers of glycoproteins in the blood plasma tend to decrease for yoga participants. The longer the participants have enrolled in the program, the greater the change in the number of detected glycoproteins for the same 4-Month of yoga Program. Glycoprotein compositions of healthy elderly people have been obtained as reference data for future related studies.

Novel Aspect

This is a comparative plasma glycoproteomics between elderly healthy yoga groups to see the efficacy of yoga programs with lectin affinity selection.

References

1. Mark D. T., Robert G. H., Jake L., Ezra A. A., *Prev Cardiol.*, 4, 165-170 (2001)
2. Wonryeon C., Kwanyoung J., Fred E. R., *J. Proteome Res.*, 9, 5960-5968 (2010)
3. Kwanyoung J., Wonryeon C., *Anal. Chem.*, 85, 7125-7132 (2013)

MP-49 / SERUM PROTEIN CLASSIFIERS OF ALS FOR DISEASE ETIOLOGY AND PROGRESSION

Szymon Filip⁽¹⁾ - *Tori Sosnowski*⁽¹⁾ - *Ina Dervishi*⁽²⁾ - *Young Ah Goo*⁽²⁾ - *Pembe Hande Ozdinler*⁽¹⁾

Northwestern University, Proteomics Center of Excellence, Chicago, United States (1) - Northwestern University, Neurology,, Chicago, United States (2)

Introduction

Amyotrophic Lateral Sclerosis (ALS) is a progressive neurodegenerative disease that affects motor neurons in the brain and the spinal cord. The progressive motor neuron degeneration eventually leads to death usually within 5 years of diagnosis [1]. ALS can be often mistaken for primary lateral sclerosis (PLS), which progresses more slowly [2]. In addition, comorbidities in ALS patients further complicates diagnosis, treatment and proper care [3, 4].

Methods

In an effort to investigate protein biomarkers for better diagnosis and disease progression, serum samples from a longitudinal cohort of 48 subjects, including ALS, re-classified from PLS to ALS, and age- and sex-matched healthy controls, were investigated. Highly abundant serum proteins were depleted and subjected to LC-MS/MS analysis. For label-free protein quantification and data analysis, MaxQuant, Perseus and R software were employed.

Results

Thirty-nine differentially expressed proteins between ALS patients (including subjects re-classified from PLS to ALS) compared to healthy controls were identified. These proteins include complement factor H, complement C2 and apolipoprotein A1, suggesting abnormality or impairment in the regulation of complement cascade and cholesterol metabolism pathways. These findings were in agreement with the previous published studies[5, 6]. Only six proteins were common between the ALS and reclassified ALS when they were compared to the healthy control. This could be due to human heterogeneity or may actually representing different molecular mechanisms affecting disease progression, which warrant further investigation. Analysis of the longitudinal patient samples revealed a number of proteins potentially

involved in disease progression. The role of candidate proteins in ALS pathophysiology and progression will be further evaluated.

Conclusions

Our study revealed the presence of a distinct set of proteins in the serum of ALS patients that can potentially be used for tracking disease progression and diagnosis. A better understanding of the pathological features characterizing ALS will shed light onto the common and unique cellular events shared among disorders of similar etiology. This will ultimately help to develop effective diagnostic tools and will potentially identify novel treatment strategies.

Novel Aspect

Identification of diagnostic protein markers for ALS disease stratification from disorders of similar etiology and for tracking ALS disease progression.

References

1. Fogarty MJ. Driven to decay: Excitability and synaptic abnormalities in amyotrophic lateral sclerosis. *Brain Res Bull.* 2018;140:318-33.
2. Singer MA, Statland JM, Wolfe GI, Barohn RJ. Primary lateral sclerosis. *Muscle Nerve.* 2007;35(3):291-302.
3. Zarei S, Carr K, Reiley L, Diaz K, Guerra O, Altamirano PF, et al. A comprehensive review of amyotrophic lateral sclerosis. *Surg Neurol Int.* 2015;6:171.
4. Fang F, Al-Chalabi A, Ronnevi LO, Turner MR, Wirdefeldt K, Kamel F, et al. Amyotrophic lateral sclerosis and cancer: a register-based study in Sweden. *Amyotroph Lateral Scler Frontotemporal Degener.* 2013; 14(5-6):362-8.
5. Lee JD, Kamaruzaman NA, Fung JNT, Taylor SM, Turner BJ, Atkin JD, et al. Dysregulation of the complement cascade in the hSOD1(G93A) transgenic mouse model of amyotrophic lateral sclerosis. *J Neuroinflamm.* 2013; 10.
6. Abdel-Khalik J, Yutuc E, Crick PJ, Gustafsson JA, Warner M, Roman G, et al. Defective cholesterol metabolism in amyotrophic lateral sclerosis. *J Lipid Res.* 2017; 58(1):267-78.

MP-50 / BLOOD PLASMA PROTEOME AND PEPTIDOME COMPARATIVE SCREENING FOR POTENTIAL BIOMARKERS OF ALZHEIMER'S DISEASE

Maria Indeykina⁽¹⁾ - ***Anna Bugrova***⁽¹⁾ - ***Natalia Zakharova***⁽¹⁾ - ***Alexey Kononikhin***⁽²⁾ - ***Yana Fedorova***⁽³⁾ - ***Alexander Brhozovskiy***⁽⁴⁾ - ***Igor Popov***⁽⁵⁾ - ***Svetlana Gavrilova***⁽³⁾ - ***Evgeny Nikolaev***⁽⁴⁾

Emanuel Institute of Biochemical Physics of the Russian Academy of Sciences, Laboratory for Mass-spectrometry of biomacromolecules, Moscow, Russian Federation (1) - Institute for Energy Problems of Chemical Physics, Laboratory for Ion and molecular physics, Moscow, Russian Federation (2) - Mental Health Research Center of the Russian Academy of Medical Sciences, Department of geriatric psychiatry, Moscow, Russian Federation (3) - Skolkovo Institute of Science and Technology, Center of Life Sciences, Skolkovo, Russian Federation (4) - Moscow Institute of Physics and Technology, Laboratory for Ion and molecular physics, Dolgoprudnij, Russian Federation (5)

Keywords: Amyloid-beta, Alzheimer's disease, blood, proteomics, peptidomics

Introduction

Alzheimer's disease is a fatal neurodegenerative disorder of the elderly, characterized by slow onset of memory loss, impaired consciousness and other brain functions, which are not detected until severe and possibly irreversible signs appear. Study of proteome changes during all stages of pathology— normal aging,

mild cognitive impairment and dementia- may allow to obtain understanding on the mechanisms of disease onset or approaches to monitor progression.

Methods

Comparative analysis of blood plasma proteins, peptides and their post-translational modifications was carried out for three groups of patients at different stages of the disorder. Proteome and peptidome fractions were extracted by a combined set of preparative procedures and analyzed by HPLC MS/MS analysis on an Agilent 1100 HPLC system coupled to a 7T FT ICR LTQ FT Ultra mass-spectrometer in DDA mode.

Results

Label-free analysis of obtained data allowed to show potential protein markers specific for different levels of cognitive disorders. Over 300 proteins were identified in total. Among the common 76 core proteins, 15 proteins were shown to reliably change in the MCI group in comparison with the controls. Of them 8 proteins were associated with various neurological diseases. A similar approach was used to find the significantly changing proteins for the AD group. A significantly less pronounced difference was observed between the AD and MCI conditions. In addition, a search for significant differences in the post-translational modifications of plasma proteins was carried out with a special attention to oxidative modifications known to be one of the early and valid tissue damage indicators, especially, for brain tissues. Significant changes in levels of PTMs such as deamidation of glutamine residues, oxidation of methionine and N-terminal residues and glutathion binding were observed.

Conclusions

The found difference in the proteomes between groups of healthy individuals and patients with MCI and a less significant variation between mild and heavy disease states allow to hope for a possibility to develop methods of diagnosis at the earlier and potentially more treatable stages of the disease progression.

Novel Aspect

The obtained protein panel can be potentially used for validation of biomarkers for screening of significant cognitive function changes and AD onset in particular.

MP-51 / A ONLINE 2D-NCFC-RP/RPLC SYSTEM FOR COMPREHENSIVE AND HIGH-THROUGHPUT PROTEOMICS ANALYSIS

Hangyeore Lee - Sang-Won Lee

Korea University, Dept. of Chemistry, Seoul, Republic of Korea

Keywords: Proteomics, UHPLC, 2DRP/RPLC, NCFC, High-throughput LC

Introduction

Proteomics aims to complete profiling of the protein contents and their modifications in cells, tissues, and biofluids and to quantitatively determine changes in their abundances. This information reveals cellular processes and signaling pathways and serves to identify candidate protein biomarkers and/or therapeutic targets. Therefore, analysis should be comprehensive and efficient.

Methods

A novel online two-dimensional reverse-phase/ reverse-phase liquid chromatography separation platform, using six switching valves based on a newly developed online non-contiguous fractionating and concatenating device (NCFC device)

Results

In bottom-up proteomics analyses of a complex proteome, this system provided significantly improved exploitation of the separation space of the two RPs, resulting in a considerable increase in the numbers of peptides identified compared to a conventional contiguous 2D-RP/RPLC method.

Conclusions

The fully automated online 2D-NCFC-RP/RPLC system bypassed a number of labor-intensive manual processes (offline fractionations, pooling, clean-up, drying/reconstitution, and autosampler injection) required with previously described offline 2D-NCFC-RP/RPLC method, offering minimal sample loss and highly reproducible 2D-RP/RPLC experiment

Novel Aspect

The developed system was the first online 2D RP/RPLC system that NCFC technology was applied. This system allowed comprehensive proteome profiling using small amount of sample without loss of sample and their information.

References

1. Yuexi W., Feng Y., Marina A. G., Yingchun W., Therese C., Tao L., Yufeng S., Matthew E. M., Daniel L., Theresa R., Ronald J. M., Richard L. K., David G. C. Richard D. S., *Proteomics*, 11, 2019-2026 (2011)
2. Hangyeore L., Dong-Gi M., Jeong Eun S., Jingi B., Hokeun K., Christophe M., Sang-Won L., *Analytical Chemistry*, 88, 23, 11734-11741(2016)

MPS-S03 - Ion spectroscopy & Photodissociation

Chairs: Nicolas C. Polfer, Simonetta Fornarini

MP-52 / STRUCTURAL ELUCIDATION OF METABOLITES USING INFRARED ION SPECTROSCOPY

Herman Bentlage⁽¹⁾ - **Jonathan Martens**⁽¹⁾ - **Giel Berden**⁽¹⁾ - **Ron Wevers**⁽²⁾ - **Jos Oomens**⁽¹⁾

FELIX, Radboud University, Nijmegen (1) - Translational Metabolic Laboratory, Radboud University Medical Center, Nijmegen (2)

Keywords: IRMPD, Metabolomics, IR-IS, FELIX, MS/MS

Introduction

The identification of molecular biomarkers is critical for diagnosis, treatment and for establishing a fundamental understanding of the pathophysiology and underlying biochemistry of inborn errors of metabolism. We show that, using the free electron laser FELIX and a modified Bruker ion trap, infrared ion spectroscopy (IR-IS) is a powerful technique for identifying metabolites. We also show that with our HPLC we can now investigate complex body fluid samples.

Methods

The experiments were performed on a modified quadrupole ion trap mass spectrometer [1] (Bruker, AmaZon Speed ETD) coupled to the infrared beam line of FELIX FEL. The FEL produces infrared radiation as 5-10 μ s macropulses at a 10 Hz repetition rate having 40-100 mJ pulse energy and a bandwidth of 0.4% of the center frequency. After one macropulse dissociation occurs. An infrared vibrational spectrum is created by relating the parent and fragment ion intensities.

Results

In a first example [2] we characterize and distinguish two isobaric biomarkers of inborn errors of

metabolism, glutaric acid for glutaric aciduria and ethylmalonic acid for short-chain acyl-CoA dehydrogenase deficiency. We show that the metabolites can be identified using an IR spectrum of a reference compound, but also in a reference-free approach by comparing the spectra against quantum-chemically calculated spectra.

In a second example [3], we show preliminary results from a project in which we attempt to identify two new biomarkers for phenylketonuria (PKU). Using an untargeted metabolomics strategy based on UPLC-qTOF mass spectrometry, several new m/z features were detected in PKU patients. However, their corresponding molecular structures could not be unequivocally assigned based only on retention time, mass and MS/MS fragmentation. Here we show our recent progress using infrared ion spectroscopy in combination with theoretical predictions and synthetic reference standards.

Conclusions

Infrared Ion Spectroscopy (IR-IS) is a powerful technique for identifying and distinguishing many kinds of molecules. Recently we have shown that the combination of liquid chromatography – mass spectrometry with IR light of FELIX creates a unique pathway to identify metabolites in complex samples. We also show that when reference standards are unavailable we can fall back on quantum-chemically calculated IR spectra for the identification of molecules.

Novel Aspect

Combining Liquid-Chromatography-Mass Spectrometry and the infrared light from FELIX we are able to identify biomarkers from complex body fluids.

References

1. Martens, J. ; Berden, G. ; Gebhardt, C. R. ; Oomens, J., Review of Scientific Instruments 2016, 87 (10), 103108.
2. Martens, J. ; Berden, G. ; Bentlage, H. ; Coene, K. L. M. ; Engelke, U. F. ; Wishart, D. ; van Scherpenzeel, M. ; Kluijtmans, L. A. J. ; Wevers, R. A. ; Oomens, J., Journal of Inherited Metabolic Disease 2018.
3. Václavík, J. ; Coene, K. L. M. ; Vrobel, I. ; Najdekr, L. ; Friedecký, D. ; Karlíková, R. ; Mádrová, L. ; Petsalo, A. ; Engelke, U. F. H. ; van Wegberg, A. ; Kluijtmans, L. A. J. ; Adam, T. ; Wevers, R. A., Journal of Inherited Metabolic Disease 2017.

MP-53 / INFRARED ION SPECTROSCOPY OF METAL-LIGAND REDOX PAIRS GENERATED VIA ELECTRON TRANSFER REDUCTION REACTIONS

Musleh Munshi - Jonathan Martens - Giel Berden - Jos Oomens

Radboud University, FELIX Laboratory, Nijmegen

Keywords: infrared ion spectroscopy, electron transfer reduction, metal-ligand complexes, redox pairs

Introduction

Transition metal organometallic complexes are known for their catalytic activation of small molecules. The catalytic activity is typically based on the ability of the metal to shuttle between different oxidation states. Here we present a novel method to spectroscopically investigate both members of a redox pair – a metal(II)-ligand complex and its charge-reduced metal(I)-ligand counterpart – in the complete isolation of the gas phase.

Methods

Infrared ion spectra were recorded using the tunable infrared free electron laser FELIX in combination with two tandem mass spectrometers, a home built Fourier-Transform Ion Cyclotron Resonance (FTICR) instrument and a Bruker AmaZon Speed ETD quadrupole ion trap (QIT) instrument [1]. The electron

transfer dissociation (ETD) option of the QIT can also be used for reducing the charge of metal(II) complexes via an electron transfer reduction (ETR) reaction with the fluoranthene radical anion [2] as electron donor.

Results

Here we demonstrate that isolated metal(I)ligand complexes can be prepared through non dissociative gas-phase one-electron reduction of the corresponding metal (II) ligand complexes. The structure of both the precursor and the charge-reduced ion, *i.e.* both members of the redox pair, were established with infrared ion spectroscopy in combination with quantum-chemical calculations. Several metal-ligand systems have been studied, including Ni-cyclam which is especially interesting for electrochemical reduction of CO₂ [3].

Conclusions

A varied set of Metal(I) ligand complexes have been generated via an electron transfer reduction reaction from the associated metal(II) complexes. The structure of these complexes in both oxidation states have been characterized via infrared ion spectroscopy using the FELIX free electron laser.

Novel Aspect

Electron transfer reduction reactions of metal (II) ligand complexes in the gas phase characterized by infrared ion spectroscopy.

References

1. Martens J., Berden G., Gebhart C. R., Oomens J., Review of Scientific Instruments, 87, 103108 (2016)
2. Martens J., Berden G., Oomens J., Analytical Chemistry, 88, 6126 (2016)
3. Munshi M. U., Craig S. M., Berden G., Martens J., De Blase A. F., Foreman J., McLuckey S. A., Oomens J., Johnson M. A., Journal of Physical Chemistry Letters, 8, 5047 (2017)

MP-55 / ION MOBILITY SPECTROMETRY AND GAS PHASE SPECTROSCOPY OF MASS AND ISOMER SELECTED PORPHYRIN ANIONS

Katrina Brendle - Patrick Weis - Manfred M. Kappes

Karlsruhe Institute of Technology, Institute of Physical Chemistry, Karlsruhe, Germany

Keywords: Porphyrins, Gas phase Spectroscopy, Ion Mobility Mass Spectrometry, Orbitrap

Introduction

Even though gas phase studies are essential for the understanding of structure-property relationships in oligomer metal porphyrins – in particular concerning their optical and chemical properties, there is only little information on isolated porphyrin ions in gas phase. Therefore we are presently investigating anionic metal porphyrin oligomers (ligand= meso-Tetra(4-sulfonatophenyl) porphine= TPPS) with different metal centers M= Cu, Zn, Pd, Mn & Fe [1,2,3].

Methods

A combination of an ion mobility drift cell for determination of ion mobilities and isomer separation, quadrupole mass filters and time-of-flight photoelectron spectrometer allows isomer and mass resolved photoelectron spectroscopy (PES) of electrosprayed analytes. Additionally we use a high resolution orbitrap mass spectrometer to identify monomer and oligomer porphyrin ions of the same m/z , which coexist in gas phase and for CID measurements.

Results

IMS measurements in combination with quantum chemical calculations have shown for [M(II)TPPS] dimers with M= Pd, Cu & Zn two different structure types [1]. They are independent of the metal center and depend strongly on charge state. Dimers in low charge states with a correspondingly larger number of

sodium counterions tend to form stacked, cofacial structures. In higher charge states, the stacks form coplanar structures. For the trivalent metal centers Mn and Fe a third “daisy-chain” like structures observed with bonds between the sulfonic acid group of one porphyrin unit and the metal center of another [2].

Photoelectron spectra upon Soret (around 420 nm) and Q band (around 550 nm) excitation of [PdTPPS]⁴⁻ reveal different electron loss channels, which were determined to correspond to direct detachment from the S1 and/or T1 state as well as tunneling emission from these two states. The first triplet state lifetime of isolated [PdTPPS]⁴⁻ is longer than 10 μs [3]. Compared to [PdTPPS]⁴⁻ the photoelectron spectra of [ZnTPPS]⁴⁻ look quite different.

Conclusions

We found different structure types depending on the total charge state and the valence of the metal center. The stacked, cofacial structure is reminiscent of H-aggregate motifs observed in solution[4].

Photoelectron spectra obtained upon excitation of the Soret and Q bands of mass selected porphyrin anions are a powerful tool to understand photoelectron detachment mechanisms. They show a strong influence of the metal center on the relaxation dynamics.

Novel Aspect

The determination of isomer and mass selected porphyrin gas phase structures is new as is the measurement of delayed electron detachment processes via isomer selective PES.

References

1. Brendle K., U. Schwarz, Jäger P., Weis P., Kappes M. M., J. Phys. Chem. A, 120, 8716 (2016).
2. Schwarz U., Vonderach M., Kappes M., Kelting R., Brendle K., Weis P., Int. J. Mass Spectrom., 339, 24 (2013)
3. Jäger P., Brendle K., Schwarz U., Himmelsbach M., Armbruster M. K., Fink K., Weis P., Kappes M. M., J. Phys. Chem. Lett., 7, 1167 (2016)
4. Ribó J. M., Crusats J., Farrera J. A., Valero M. L., J. Chem. Soc., Chem. Commun., 681 (1994)

MP-56 / REMPI AND MATI SPECTROSCOPY OF METHYL AND HALO PYRIDINES

Niklas Helle - Jürgen Grotemeyer

Institute of Physical Chemistry, CAU Kiel, Kiel

Keywords: REMPI, MATI, ionization energies, vibronic structure, pyridine

Introduction

Pyridine has been extensively investigated by means of REMPI and MATI spectroscopy [1,2]. In contrast, only little information can be found about the vibronic structure of the methyl and halo pyridines in their first electronically excited and ionic ground state. Hence, we focused on analyzing the electronic and geometrical structure of the first excited state S1 ($1\pi^* \leftarrow n$) and ionic ground state D0 of the picolines [3], 3-fluoro and 3-chloropyridine.

Methods

The molecules were investigated by means of REMPI and MATI spectroscopy. The measurements were done with two dye lasers and a home-built ToFMS.

In REMPI mode, the spectra were obtained by scanning the excitation laser. In MATI mode, the excitation laser was kept resonant with an eigenstate and the ionization laser was scanned. The immediately formed ions were separated from the Rydberg neutrals by a spoiling field. Afterwards, the Rydberg molecules were field ionized.

Results

The excitation energy of the picolines, 3-fluoro and 3-chloropyridine have been determined with a very high accuracy. The influence of the position of the methyl group and the halogen atoms has been identified. All derivatives experience an in plane geometry distortion along the eigenvector of the vibration 6a upon S1 excitation. Except for 3-chloropyridine, this vibration has the highest contribution to the electronic transition. Furthermore, an out of plane geometry distortion along the eigenvector of low frequency vibrations, which have a varying intensity among the derivatives, has been observed.

The adiabatic ionization energy of the picolines, 3-fluoro and 3-chloropyridine has been determined by means of MATI spectroscopy. The picolines and 3-fluoro pyridine experience a geometry distortion along the eigenvector of the vibration 8b upon ionization via the S1 state. In contrast, the corresponding MATI spectrum of 3-chloro pyridine is dominated by low frequency vibrations.

Conclusions

The ionization energy of the pyridine derivatives were measured with the highest accuracy reported so far. In addition, the ionic structure of the molecules has been determined.

The MATI spectra via the active vibrations in the S1 state offered the possibility to verify the assignment in the S1 REMPI spectra. This was especially true for low frequency vibrations, which are symmetry forbidden and were not taken into account for the simulations of the vibrational spectra.

Novel Aspect

The adiabatic ionization energies of the picolines, 3-fluoro and 3-chloro pyridine were measured with a very high accuracy and the ionic geometries were determined.

References

1. Riese M., Altug Z., Grotemeyer J., Phys. Chem. Chem. Phys., 8, 4441-4448 (2006).
2. Lee Y. R., Kang D. W., Kim H. L., Kwon C. H., J. Chem. Phys., 141, 174303 (2014).
3. Helle N., Krüger S., Grotemeyer J., Curr Phys Chem, 8, in press (2018).

MP-57 / ACTION SPECTROSCOPY WITHOUT A TRAP: ON-THE-FLY SPECTROSCOPY OF IR-797 AND IR-806 CYANINE DYES

Marius Herve - Richard Bredy - Gabriel Karras - Bruno Concina - Abdul Rahman Allouche - Franck Lepine - Isabelle Compagnon

Institut Lumière Matière (ILM), Université Lyon 1, Lyon

Keywords: Cyanine dyes (IR-797, IR-806); Mass Spectrometry; NIR Action Spectroscopy; Triple Quadrupole; Geometry sensitive absorption shift

Introduction

Ion spectroscopy has become a powerful method in the modern MS toolkit. Numerous experimental setups make use of ion traps to irradiate molecules, thus accessing absorption of ions in the gas phase. We used here a triple quadrupole instrument to measure photodissociation action spectra of IR-797, a heptamethine cyanine dye, and one of its doubly sulfonated derivative IR-806 [1]. This revealed absorption shifts depending on the sulfonate presence and counter cation.

Methods

Gas phase cyanines were produced using electrospray ionization, and transferred into a triple quadrupole mass spectrometer (Xevo TQ-S micro, Waters). Ions selected in charge and counter cation state were irradiated after the first mass selecting quadrupole with tunable femtosecond laser pulses from 710 nm to 900 nm (80 MHz, 140 fs, 25 nJ/pulse, Coherent Chameleon Ultra II). Photofragments were then analyzed by the last quadrupole.

Results

Photofragmentation yield was recorded as a function of the laser wavelength resulting in NIR photodissociation action spectra. Maximum photodissociation of $[\text{IR-797}]^+$ is observed at 760 nm, which is blue-shifted compared to the maximum absorption in the aqueous phase, and corresponds to the S0-S1 transition. The addition of sulfonated substituents on the nitrogen atoms leads to an increase of the maximum photodissociation yield to 780 nm for negative $[\text{IR-806}]^-$ ions and doubly protonated $[\text{IR-806,2H}]^+$ ions. This results from the modification of the geometry of the cyanine chromophore with sulfonate groups in the vicinity. We demonstrated this effect by changing the sulfonate counter cation. Changing the counter cation from H^+ to Na^+ further modifies the position of the sulfonate substituents around the chromophore, yielding a supplementary absorption shift to 790 nm.

Power dependence also revealed that two photons are necessary to fragment, implying sequential absorption to the first excited state S1 followed by intramolecular relaxation.

Conclusions

Photodissociation action spectroscopy of gas phase IR-797 and IR-806 cyanine ions reported shifts in the absorption to the S1 excited state. These shifts constitute a probe of the change in the geometry of the molecule, where the absorption of the chromophore reveals its close environment composed of sulfonate with or without counter cation [2]. This can have consequences for aqueous phase absorption in the presence of ion pairs [3].

Novel Aspect

We performed counter cation-dependent action spectroscopy on cyanine dyes using a triple quadrupole instrument coupled to a tunable high-repetition rate femtosecond laser.

References

1. Nairat M., Webb M., Esch, M., Lozovoy V., Levine B., Dantus M., Phys. Chem. Chem. Phys., 19, 14085 (2017).
2. Stockett M., Boesen M., Houmøller J., Nielsen S., Angew. Chem. Int. Ed., 56, 3490-3495 (2017)
3. Tatikolov A., Shvedova L., Derevyanko N., Ishchenko A., Kuzmin V., Chem. Phys. Lett., 190, 291-297 (1992)

MP-58 / IR TAGGING PRE-DISSOCIATION SPECTROSCOPY TO REVEAL ISOMERIC COMPOSITION

Pavol Jusko⁽¹⁾ - *Aude Simon*⁽²⁾ - *Gabi Wenzel*⁽¹⁾ - *Sandra Brünken*⁽³⁾ - *Stephan Schlemmer*⁽⁴⁾ - *Christine Joblin*⁽¹⁾

Institut de Recherche en Astrophysique et Planétologie (IRAP), Université de Toulouse, CNRS, CNES, Toulouse (1) - LCPO, Université de Toulouse, CNRS, Toulouse (2) - FELIX, Radboud University, Nijmegen (3) - 1. Physikalisches Institut, Universität zu Köln, Cologne (4)

Keywords: Cryogenic ion trap, tagging spectroscopy, isomer analysis, tropylium, benzylium

Introduction

The recent introduction of cryogenic ion traps reaching temperatures of about 4 K opens new possibilities in weak bond ion-neutral attachment, notably the possibility of He/Ne attachment to the ions of interest [1]. The infrared pre-dissociation of these tagged ions can be used in ro-vibrational action spectroscopy. This can be used to infer the isomer composition of the ions.

Methods

The FELion beamline at FELIX free electron laser facility at Radboud University (Nijmegen, The Netherlands), is equipped with a 22 pole radio-frequency ion trap at 4 K [2]. The ions, produced by electron bombardment are mass selected prior to injection into the trap, where they are tagged with He or Ne.

Consecutively tunable FEL light is used to irradiate the tagged ions and obtain a one photon action spectrum (IR pre-dissociation).

Results

The tropylium (7 carbon-hydrogen aromatic cycle) and benzylium ions (6 carbon-hydrogen aromatic cycle with one hydrogen substituted by methylene CH_2 group) are the two lowest-energy isomers of C_7H_7^+ . Tropylium and benzylium like structures could also be found in larger polycyclic aromatic hydrocarbon (PAH) species, *e.g.*, in a fragment of 1-methylpyrene $\text{C}_{17}\text{H}_{11}^+$ (denoted as PyrC7^+ and PyrCH_2^+) as predicted by Rapacioli et al. [3]. In this study, we will demonstrate how the IR pre-dissociation spectroscopy of cold Ne-tagged ions can reveal the isomeric mixture in situations where the IRMPD spectroscopy fails. We applied molecular dynamics simulations (SCC-DFTB level of theory) in order to analyse the difference between the multiphoton and single photon spectra. We also compare our results with isomer mixture ratios that can be obtained from the variation in reactivity of the two isomers.

Conclusions

We show, that IR pre-dissociation of the weakly bound Ne cluster can be used as a complementary technique to mass spectrometry, providing information on isomeric composition. We identified the most abundant isomer of cationic 1-methylpyrene fragment, produced by low energy electron bombardment of the neutral 1-methylpyrene, to be PyrCH_2^+ isomer [4].

Novel Aspect

The single photon dissociation scheme has a great advantage compared to multi-photon ones, in which consecutive heating can lead to isomerisation.

References

1. Chakrabarty S., Holz M., Campbell E. K., Banerjee A., Gerlich D., Maier J. P., J. Phys. Chem. Lett., 4, 4051 (2013).
2. Jusko P., Brünken S., Asvny O., Thorwirth S., Stoffels A., Redlich B., van der Meer L., Oomens J., Schlemmer S., In prep.
3. Rapacioli M., Simon A., Marshall C. C. M., Cuny J., Kokkin D., Spiegelman F., Joblin C., J. Phys. Chem. A, 119, 12845 (2015).
4. Jusko P., Simon A., Wenzel G., Brünken S., Schlemmer S., Joblin C., Chem. Phys. Lett., 698, 206 (2018).

MP-59 / VALIDATION OF AN IRMPD EXPERIMENTAL SETUP AND THEORETICAL METHODS FOR IRMPD ASSIGNMENT

André F. Rodrigues-Oliveira - Francisco W. M. Ribeiro - Gustavo Cervi - Thiago C. Correra

Institute of Chemistry, University of São Paulo, São Paulo

Keywords: IRMPD; theoretical calculations; intramolecular interaction; ion spectroscopy; bioactive molecules

Introduction

Ion spectroscopy are usually compared to theoretical spectra to determine conformation of gaseous ions. DFT methods, mainly B3LYP, are the most used methods to model IR spectra, but some of them are known to provide inaccurate description of highly anharmonic oscillators. [1,2] We aim to assess alternative methods for describing anharmonic oscillators by comparison with experimental IRMPD spectra.

Methods

The benchmarking and validation of our experimental setup was performed with experimental spectra of protonated biotin and tryptophan. Photodissociation was achieved by an OPO/A system in the 2750 – 3800

cm⁻¹ range. Calculations were performed with B3LYP, CAM-B3LYP, M06-2X and MP2 and the basis sets of Pople and Dunning type. Anharmonicity and D3 dispersion corrections were also evaluated.

Results

M06-2X/6-31+G** was the most cost effective level of theory. The increase in the job time, compared to B3LYP, was not significant, while the error (RMSD) highly decreases. Extended basis sets were not as much efficient, when associated to M06-2X. CAM-B3LYP was the better method to describe intramolecular interactions and provided good results with extended Dunning type basis sets. Anharmonic calculations were carried out and, but they showed to be extremely expensive when compared to other approaches. The dispersion correction D3, improved the description of O-H vibrations involved in H-Bond, but worsened the description of N-H stretches. None of the correction methods lowered RMSD as the use of well established scaling factors (sf). The M06-2X/6-31+G** and CAM-B3LYP/aug-cc-pVDZ level of theory was also tested against the fingerprint region, showing a good agreement with experimental data.

Conclusions

Our results suggest that M06-2X/6-31+G(d,p) could be the level of theory to be used in predicting IRMPD spectra. On the other hand, CAM-B3LYP/aug-cc-pVDZ was more efficient in describing long-range interactions, despite requiring more time than the previous level of theory. Between the correction approaches, the use of scaling factors remains as the easiest and more efficient strategy.

Novel Aspect

These are first IRMPD done in Brazil. We show the most used levels of theory to predict IRMPD spectra. We suggest the use of M06-2X/6-31+G**.

References

1. M. G. Medvedev, I. S. Bushmarinov, J. Sun, J. P. Perdew and K. A. Lyssenko, *Science*, 355, 49–52, 2017
2. A. J. Cohen, P. Mori-Sanchez and W. Yang, *Science*, 321, 792–794, 2008.

MP-60 / INFRARED MULTIPLE PHOTON DISSOCIATION OF COPPER FORMATE CLUSTERS

Tobias Pascher⁽¹⁾ - ***Milan Oncak***⁽¹⁾ - ***Christian van der Linde***⁽¹⁾ - ***Josef Lengyel***⁽²⁾ - ***Martin K. Beyer***⁽¹⁾

Universität Innsbruck, Institut für Ionenphysik und Angewandte Physik, Innsbruck (1) - Technische Universität München, Physikalische Chemie, München (3)

Keywords: FTMS, Catalysis, Copper, Formate, IRMPD, Spectroscopy

Introduction

Copper and copper oxides are widely investigated for their catalytic properties [1,2]. It was shown that the calcination of copper salts can lead to highly active catalysts [3] and that copper hydride-based catalysts show a distinct reactivity [4]. Here, elementary steps in the calcination of copper formate are investigated in the gas phase by means of infrared multiple photon dissociation of $\text{Cu}_n(\text{HCOO})_{2n+1}^-$, $n \leq 4$.

Methods

$\text{Cu}_n(\text{HCO}_2)_m^-$ ($n \leq 4$; $m \leq 9$) clusters are prepared using electrospray ionization, transferred into the cell of a Fourier transformation ion cyclotron resonance mass spectrometer and irradiated in the wavelength range of 2500-4475 nm by a tuneable optical parametric oscillator. The decomposition of the ions is investigated by measuring dissociation kinetics at the absorption bands. For theoretical comparison, quantum chemical calculations on the DFT level are used.

Results

Using action spectroscopy, three intense absorption bands were found in the wavelength range of 3400-

3900 nm that shift depending on the size of the cluster and the oxidation state of the copper ions. Using deuterated isotopologues, the character of the vibrations involving C-H/D motion was confirmed and absorptions of C-O combination bands/overtone were identified.

After activating the C-H bonds of the formate ligands with the tuneable laser system, a fragmentation cascade has been observed with heavy secondary fragmentation. Clusters $\text{Cu}_n(\text{HCOO})_{2n+1}^-$, $n = 2-4$, dissociate into smaller $\text{Cu}_n(\text{HCOO})_{2n+1}^-$ species with $n = 1-3$. For $\text{Cu}_2(\text{HCOO})_5^-$, H atom transfer takes place between two formate ligands, resulting in the elimination of CO_2 and HCOOH together with a reduction of the copper centers to +I. For smaller clusters ($n = 1,2$; $m = 2,3$), CO_2 elimination involves the transfer of a hydrogen atom to copper, leading to copper hydrides. Collisional activation leads to the same fragmentation pathway.

Conclusions

Three intense absorption bands of C-H stretching vibrations were found in the wavelength range of 3400-3900 nm that shift depending on the size of the cluster and the oxidation state of the copper ions. The dissociation of copper formate ions was found to lead towards copper hydrides through a reaction cascade. The results were compared to quantum chemical calculations on the DFT level.

Novel Aspect

Copper formate clusters are activated with infrared laser light to simulate calcination with defined heating. CO_2 elimination leads to copper hydride species.

References

1. Elwell C. E. et al., Chem. Rev., 117, 2059–2107 (2017).
2. Tsybizova A., Roithová J., Mass Spectrom. Rev., 35, 85-110 (2016).
3. Prasad R., Singh P., Catal. Sci. Technol., 12, 3326(2013).
4. Jordan A. J., Lalic G., Sadighi J. P., Chem. Rev., 116, 8318-8372(2016).

MP-61 / DEVELOPMENT OF A 3D ION TRAP MASS SPECTROMETER FOR IRMPD SPECTROSCOPY

Tatiana C. Penna - Gustavo Cervi - André F. Rodrigues-Oliveira - Jair Menegon - Thiago C. Correra

Instituto de Química, Universidade de São Paulo, São Paulo

Keywords: IRMPD spectroscopy, IRLAPS, CO_2 laser, photofragmentation

Introduction

Ion spectroscopy experiments are excellent tools to obtain complementary information in mass spectrometry. We describe the development of the 3D ion trap used for IRMPD spectroscopy and IRLAPS experiments, comparing results already present in literature and different ionization sources. This equipment is the first to perform IRMPD experiments in Brazil and is open to use by external users.

Methods

Experiments were made in a modified Bruker AmaZon SL3D ion trap mass spectrometer. IRMPD experiments used an Laser Vision OPO/A system pumped by Nd:YAG 1064 nm laser to obtain tunable radiation from ca. 2300 – 4000 cm^{-1} . Crystal optimization were performed when necessary to improve laser power in different wavenumber ranges.

Results

A LabView interface was developed to annotate the actual laser wavelength to the raw MS data, allowing the direct calculation of the photodissociation efficiency. This interface completely controls the MS events and allows for multiple and variable IRMPD steps by customizing CID events.

IRMPD spectra of species previously reported in literature were obtained to verify the instrument performance using ESI-MS. The spectra for biotin, tryptophan, phosphotyrosine were similar to those previously reported which indicates the good performance of our system. [1-3]

IRMPD of the $\text{Ag}^+(\text{CH}_3\text{CN})_2$ species were obtained at the nitrile stretch region as dissociation at 2350 cm^{-1} was registered, which is unprecedented given the low output power of the OPO/A at this wavelength.

IRMPD spectra of a series of ions were obtained by using other in house built ion sources like, nanospray, EASI and DESI ion sources.

Conclusions

We have demonstrated that IRMPD spectra obtained in our new instrument are consistent with previous results of the literature. It was also possible to obtain the vibrational spectra from different ionization sources. The optimization of the OPO crystals positions improved the output power of the pulse energy in different ranges of wavenumber and allowed the IRMPD of nitriles at 2350 cm^{-1} .

Novel Aspect

The recently assembled instrument is the first to perform IRMPD experiments in Brazil. IRMPD were obtained from different ion sources and at 2350 cm^{-1} .

References

1. Fraschetti, C., Filippi, A., Guarcini, L., Steinmetz, V., Speranza, M., J. Phys Chem. B, 119, 6198 (2015).
2. Mino, W. K., Gulyuz, K., Wang, D., Stedwell, C. N., Polfer, N. C., J. Phys. Chem. Lett., 2, 299 (2011).
3. Scuderi, D., Bakker, J. M., Durand, S., Maitre, P., Sharma, A., Martens, J. K., Nicol, E., Int. J. Mass Spec., 308, 338 (2011).

MP-62 / PHOTOISOMERIZATION ACTION SPECTROSCOPY OF A RUTHENIUM SULFOXIDE COMPLEX

Michael Scholz⁽¹⁾ - ***James Bull***⁽¹⁾ - ***Eduardo Carrascosa***⁽¹⁾ - ***Gilbert Kosgei***⁽²⁾ - ***Jeffrey Rack***⁽²⁾ - ***Evan Bieske***⁽¹⁾

University of Melbourne, School of Chemistry, Parkville, Australia (1) - University of New Mexico, Department of Chemistry and Chemical Biology, Albuquerque, United States (2)

Keywords: ion mobility spectrometry; ruthenium; photoisomerization; laser spectroscopy; dynamics

Introduction

Understanding the excited-state dynamics of photoactive molecules is essential for the design of new generation photoresponsive materials. One promising family of switching units is based on the reversible linkage isomerism of sulfoxide ligands coordinated to a ruthenium(II) polypyridyl core. In these, the sulfoxide ligands convert between S- and O-bound configurations with quantum yields of up to 79% after metal-to-ligand charge transfer (MLCT) excitation.

Methods

Photoisomerization action (PISA) spectroscopy was performed using a custom tandem ion mobility mass spectrometer coupled with a tunable laser. In an experiment, electrosprayed isomers were mobility separated in a first drift region, selected with an ion gate, and irradiated with tunable light. Photo-isomers were mobility separated in a second drift region and mass selected with a quadrupole. Resulting PISA spectra plot photo-isomer yield against wavelength.

Results

The $[\text{Ru}(\text{bpy})_2(\text{bpSO})]^{2+}$ complex (bpy is 2,2'-bipyridine, bpSO is a bidentate bis-sulfoxide ligand) displays stepwise isomerization in solution from $S,S \rightarrow S,O \rightarrow O,O$ with an overall quantum yield of 42% at 355 nm. [1] In the gas-phase, we observed direct S,S to O,O photoisomerization with minor formation of the S,O

isomer. [2] Supporting electronic structure calculations suggested that both $\pi\pi^*$ and $\pi\pi^*$ excitation resulted in isomerization. The difference in photoisomerization yields between the previous measurements in solution and our gas-phase results was attributed to solvent quenching of vibrational energy before activated molecules could traverse excited-state isomerization barriers.

Conclusions

Ion mobility mass spectrometry is a powerful technique for separating linkage isomers of ruthenium sulfoxide coordination complexes. Different photoisomerization yields were observed in the gas phase and in solution, indicating that solvation strongly perturbs isomerization process. Future studies will aim to elucidate further details of the photochemical dynamics in these systems.

Novel Aspect

Ion mobility mass spectrometry coupled with laser spectroscopy characterizes the first example of linkage photoisomerization in a gas-phase coordination complex.

References

1. King, A. W., Jin, Y., Engle, J. T., Ziegler, C. J., Rack, J. J., *Inorg. Chem.*, 52 (4), 2086–2093 (2013).
2. Scholz, M. S., Bull, J. N., Carrascosa, E., Adamson, B. D., Kosgei, G. K., Rack, J. J., Bieske, E. J., *Inorg. Chem.*, 57 (9), 5701–5706 (2018).

MP-63 / PHOTOCHEMISTRY OF $\text{Mg}^+(\text{H}_2\text{O})_n$, $n=1-5$: EXPERIMENT AND THEORY

Thomas Taxer - Milan Ončák - Christian van der Linde - Erik Barwa - Martin K. Beyer

Universität Innsbruck, Institut für Ionenphysik und Angewandte Physik, Innsbruck

Keywords: Photochemical hydrogen evolution; photodissociation

Introduction

There has been a long history of studies [1-4] on hydrated metal ions in the gas phase over the last decades. The study of such cluster systems is important to understand the transition of various properties from single molecules to the bulk. In the ground state these systems are already well understood but few studies are available on their excited state reactivity. Here, we investigate the photochemistry and spectroscopy of hydrated magnesium ions $\text{Mg}^+(\text{H}_2\text{O})_n$, ($n=1-5$).

Methods

The experiments were performed on an FT-ICR-MS, using a laser vaporization ion source. Clusters were generated from a target consisting of isotopically enriched ^{24}Mg and supersonic expansion into high vacuum, along with a gas pulse of helium seeded with water vapor. A tunable, pulsed laser system was used to induce photodissociation. Dissociation mass spectra were recorded in a wavelength range from 240 nm up to 2,000 nm and compared with ab initio calculations.

Results

Photodissociation spectra for $\text{Mg}^+(\text{H}_2\text{O})_n$, $n = 1-5$ were measured, at 130 ± 20 K and compared to experiments by Misaizu et al. [2]. The detected fragment ions can be grouped into two different channels. The $\text{Mg}^+(\text{H}_2\text{O})_m$ ($m < n$) species are produced by the evaporation of one or more water molecules. The other channel, contributing most to the observed fragmentation, is the hydrogen dissociation reaction producing magnesium hydroxide $\text{MgOH}^+(\text{H}_2\text{O})_m$. The additional loss of one or more water molecules was also observed in this dissociation channel. For the $\text{Mg}^+(\text{H}_2\text{O})_{4,5}$ clusters, the influence of black-body infrared radiative dissociation was documented. Detailed photofragment branching ratios were evaluated. The main bands in the experimental spectra were assigned to the $3s \rightarrow 3p$ transition. Additional bands were observed at higher photon energies, corresponding to the $3s \rightarrow 3d/4s$ transitions. Resolved vibrational

structure was obtained for $\text{Mg}^+(\text{H}_2\text{O})$. For this system the possibility of a two photon process in the dissociation was investigated theoretically.

Conclusions

Photodissociation spectra were measured for $\text{Mg}^+(\text{H}_2\text{O})_{1-5}$. The results are in good agreement with theoretical predictions and earlier experiments[2]. The $3s \rightarrow 3d/4s$ dissociation bands were observed for $\text{Mg}^+(\text{H}_2\text{O})_{3-5}$. For $\text{Mg}^+(\text{H}_2\text{O})_{1-3}$, only a single isomer was present, but several isomers contribute to the spectrum of $\text{Mg}^+(\text{H}_2\text{O})_{4,5}$. For all investigated cluster sizes, hydrogen dissociation producing $\text{MgOH}^+(\text{H}_2\text{O})_m$ was the main observed dissociation channel.

Novel Aspect

Photodissociation spectra for $\text{Mg}^+(\text{H}_2\text{O})_{1-5}$ were measured with high resolution, far into the UV, in a temperature controlled environment.

References

1. Yeh C. S., Willey K. F., Robbins D. L., Duncan M. A., Int. J. Mass Spectrom, 131, 307–317 (1994).
2. Misaizu F., Sanekata M., Fuke K., Iwata, S., J. Chem. Phys. 100, 1161–1170 (1994).
3. Berg C., Achatz U., Beyer M., Joos S., Albert G., Schindler T., Niedner-Schatteburg G., Bondybey V. E., Int. J. Mass Spectrom. 167, 723-734 (1997).
4. Berg, C. ; Beyer, M. ; Achatz, U. ; Joos, S. ; Niedner-Schatteburg, G. ; Bondybey, V. E, Chem. Phys. 239, 379–392 (1998).

MP-64 / SHEDDING LIGHT ON STRUCTURE: INITIAL PROTEIN UNFOLDING EVENTS REVEALED BY 213 NM UVPD

Alina Theisen⁽¹⁾ - Rachelle Black⁽¹⁾ - Davide Corinti⁽²⁾ - Jeffrey Brown⁽³⁾ - Bruno Bellina⁽¹⁾ - Perdita Barran⁽¹⁾

University of Manchester, Manchester Institute of Biotechnology and Photon Science Institute, Manchester (1) - Universita di Roma "La Sapienza", Dipartimento di Chimica e Tecnologie del Farmaco, Rome (2) - Waters, R&D, Wilmslow (3)

Keywords: 213 nm UVPD, Native MS, Ion mobility, Protein unfolding

Introduction

Protein unfolding in its initial stages may be reflective of the stability of the entire fold and reveal local sites that may be perturbed. Ion mobility mass spectrometry (IM-MS) provides insight into global structure while UV photodissociation (UVPD) has been shown to be sensitive to secondary structural elements [1,2,3]. In this work we investigate the initial unfolding of three proteins using 213 nm UVPD combined with IM-MS.

Methods

Modifications were made to a Synapt G2-S [4] to facilitate injection of a laser beam into the trap cell region. Ubiquitin, cytochrome c and myoglobin were sprayed from ammonium acetate by nanoESI; UVPD spectra using 213 nm photons and arrival times of precursor and photofragment ions were recorded at different source conditions by varying the cone voltage affecting the protein conformation prior to UVPD.

Results

Native structures and how they alter upon gentle in source activation was assessed. Unfolding from compact to extended conformer was induced in the lowest charge state of ubiquitin and UVPD was performed as a function of activation; the fragmentation yield was higher for the extended conformation. Cleavage in the compact structure was enhanced between residues 12-18 compared to residues 33-47 for the extended structure.

For cytochrome c UVPD was carried out on a compact, an intermediate and an extended conformation. Again the fragmentation yield increased as the protein increased in CCS with distinct UVPD spectra for each conformation; cleavage at residues Phe10, Thr19 and Val20 was only found in activating conditions. Compact conformations of myoglobin yielded few fragments while a drastic increase was observed after a shift to intermediate and extended structures. The cleavage sites with the highest yield were similar for both native charge states, however fragmentation in the GH loop, H helix and D helix was only found for 9+.

Conclusions

UVPD-IM-MS discriminates between compact structures and more extended forms. We present differences in terms of fragmentation yield as well as the intensity distribution within the fragmentation pattern. Cleavage sites differ between conformational ensembles revealing areas of local structural changes as a result of progressive in-source activation as well as global changes in terms of ion mobility profiles.

Novel Aspect

213 nm UVPD-IM-MS is used to identify and locate structural changes upon protein unfolding

References

1. Warnke S., von Helden G., Pagel K., *Proteomics*, 15, 2804-2812 (2010)
2. Cammarata M. B., Brodbelt J. S., *Chem. Sci.*, 6, 1324-1333 (2015)
3. Theisen A., Yan B., Brown J. M., Morris M., Bellina B., Barran P. E., *Anal. Chem.*, 88, 9963-9971 (2016)
4. Mistarz U. H., Bellina B., Jensen P. F., Brown J. M., Barran P. E., *Anal. Chem.*, 90, 3079-3082 (2018)

MP-65 / IRMPD SPECTRA OF PROTONATED HYDROXYBENZALDEHYDES : THE EVIDENCE OF THE TORSIONAL BARRIERS IN CARBOXONIUM IONS

Barbara Chiavarino⁽¹⁾ - ***Maria Elisa Crestoni***⁽¹⁾ - ***Otto Dopfer***⁽²⁾ - ***Philippe Maitre***⁽³⁾ - ***Simonetta Fornarini***⁽¹⁾

Università di Roma "La Sapienza", Dipartimento di Chimica e Tecnologie del Farmaco, Roma, Italy (1) - Technische Universität Berlin, Institut für Optik und Atomare Physik, Berlin, Germany (2) - Université Paris Sud,, Laboratoire de Chimie Physique UMR8000 CNRS, Orsay, France (3)

Keywords: Protonated aromatic compounds, carboxonium ions , IRMPD spectroscopy, DFT calculations,

Introduction

The basicity of aromatic aldehydes was widely studied. Protonation of benzaldehyde is directed to oxygen yielding a carboxonium ion with bivalent behavior as oxonium and carbenium ion. [1] In substituted benzaldehydes, the prevalence of the carbenium or oxonium character depends on the substituents and on their position. The aim of this work is to characterize the structure of protonated hydroxybenzaldehydes, identified as key inhibitors in biomass hydrolysates.

Methods

Protonated ortho, meta, and para-hydroxybenzaldehydes (IBZH⁺) were generated in the gas phase by electrospray ionization of hydro-alcoholic solutions, (10⁻⁵ M) of the select isomer, containing 2% acetic acid. Their structure was probed by combining IRMPD spectroscopy in the X-H spectral region with DFT calculations at the B3LYP/6-311++G(2df2p) level. As further inquiry, IRMPD spectra of both p-IBZH⁺ and its fragment ion were recorded in the fingerprint region.

Results

Similar to protonated aliphatic aldehydes, protonation at the formyl oxygen atom of benzaldehydes leading to the formation of carboxonium ions yields two distinct isomers, depending on the relative orientation of the proton that can stay either in cis or trans with respect to hydrogen on the adjacent carbon [2,3]. In this

context, the 3200-3700 cm^{-1} spectral region results ideal to distinguish the characteristic O–H stretch frequency of either cis and trans carboxonium isomers ions due to the different O–H bond properties of the two isomers. Comparison between IRMPD spectra and DFT calculations suggests that in the case of both p-IBZH⁺ and m-IBZH⁺ only cis conformers are present in the ionic population analyzed. In the case of o-IBZH⁺, IRMPD spectroscopy points to a mixture of one trans and two cis conformers. The value of energy barrier for the cis-trans isomerization calculated for each IBZH⁺ appears to be a measure of the degree of π -electron delocalization in the considered system.

Conclusions

The present study shows that protonation of hydroxybenzaldehydes is directed only to the formyl oxygen atom to generate a mixture of carboxonium ions, with exclusively cis conformation. Only in the case of protonated o-IBZH⁺, the cis carboxonium ion isomer may overcome the high rotational barrier, and gives partially the more stable trans conformer. Moreover the IRMPD spectra of a p-IBZH⁺ fragment formed by loss of CO corresponds to para-protonated phenol.

Novel Aspect

This work is the first gas-phase experimental study of IBZH⁺ and reports for the first time an IR spectrum of gaseous protonated phenol in the fingerprint region.

References

1. Olah G. A., Rasul, G., York C., Prakash G. K. S., J. Am. Chem. Soc., 117,11211–11214 (1995)
2. Chiavarino B., Crestoni M. E., Fornarini S., Dopfer O., Lemaire J., Maître P. J. Phys. Chem. A, 110, 9352-9360(2006)
3. Chakraborty S., Patzer A., Dopfer O. The Journal of Chemical Physics, 133, 044307 (2010)

MP-66 / VALENCE AND CORE ELECTRON INDUCED PHOTOIONIZATION AND DISSOCIATION DYNAMICS OF THE THIAZOLE MOLECULE

Alexandre Lago⁽¹⁾ - *Robin Januario*⁽¹⁾ - *Reinaldo Cavasso Filho*⁽¹⁾ - *Marc Simon*⁽²⁾ - *Juan Dávalos*⁽³⁾

Federal University of ABC (UFABC), Center for Natural Sciences and Humanities (CCNH), Santo Andre - SP (1) - UPMC-Univ Paris 06, CNRS, Sorbonne Universités, UPMC-Univ Paris 06 / Laboratoire de Chimie Physique-Matiere et Rayonnement, Paris (2) - Consejo Superior de Investigaciones Científicas (CSIC), CSIC / Instituto de Química Física Rocasolano, Madrid (3)

Keywords: Photoionization, Ionic dissociation, Synchrotron radiation, PEPICO, Thiazole.

Introduction

The thiazole ($\text{C}_3\text{H}_3\text{NS}$) molecule belongs to the group of organic compounds called heterocycles, which are present in the formulation of important complex chemical entities, such as carbohydrates, proteins, among others, and stand out for their diverse properties and applications important to several fields. The main goal in this project is to provide reliable fundamental spectroscopic and spectrometric data for this group of molecules in gas phase by using synchrotron radiation, photoelectron photoion coincidence and time of flight mass spectrometry techniques.

Methods

TOF-MS and photoelectron photoion coincidence (PEPICO) techniques have been employed. VUV and soft X-ray synchrotron radiation sources were employed as ionizing agents. Total and partial ion yields (TIY,PIY) and PEPICO spectra were recorded as a function of the incident photon energy, covering from the valence region up to the core electronic shells of the molecule. The resulting photoionization ion products and the possible fragmentation pathways leading to those species are presented and discussed.

Results

For the valence photodissociation studies of thiazole, the PEPICO mass spectra have been recorded covering the range from 10 eV to 21 eV. The high photo-stability of this molecule and its molecular parent ion has been confirmed from our results. It has been mainly attributed to the strong bonding resulted from its closed ring molecular structure. The dominant mechanism corresponds to the consecutive dissociations of the C=N and C=C bonds from the thiazole molecular ion resulting in the fragment $C_2H_2S^+$ ($m/z = 58$). The ion yield spectra were recorded by scanning the incident photon energies in the vicinity of the C1s and S1s core edges showing the main spectral features resulting from core electronic excitation transitions from the molecular ground state leading to the core excited and to the corresponding ionic states. In addition, the corresponding mass spectra based on the PEPICO technique recorded at selected photon energies near the core ionization edges provided relevant information leading to the identification of the dissociative photoionization ion products and their corresponding fragmentation pathways.

Conclusions

Time of flight mass spectrometry (TOF-MS) used in conjunction with photoelectron photoion coincidence (PEPICO), total and partial ion yield (TIY and PIY) techniques have been proved powerful tools for the analysis of the electronic structure and photodissociation dynamics of the thiazole molecule in gas phase, which resulted in very interesting chemical and physical insights.

References

1. A. F. Lago, et al., *J. Mass Spectrom.* 2014, 49, 1163–1170.
2. A. F. Lago, et al., *J. Mass Spectrom.* 2017, 52, 657–663.

Acknowledgements: CNPq, FAPESP, CAPES.

MP-67 / RELEVANCE OF THE GAS-PHASE STUDIES TO BIOLOGY: COLD ION SPECTROSCOPY OF THE OPIOID PEPTIDE DRUGS

Aleksandr Y. Pereverzev⁽¹⁾ - **Vladimir N. Kopysov**⁽¹⁾ - **Istvan Szabo**⁽²⁾ - **Edina Rosta**⁽²⁾ - **Oleg V. Boyarkin**⁽¹⁾

EPFL, LCPM, Lausanne, Switzerland (1) - Kings College, Department of Chemistry, London, United Kingdom (2)

Keywords: UV& IR spectroscopy, IR-IR-UV hole-burning, enkephalin, photodissociation

Cold Ion Spectroscopy (CIS) provides accurate spectroscopic fingerprints for stringent validation of the 3D structures of small biomolecules in the gas phase [1-3]. A fundamental question that still remains open is to what extent the gas-phase structures of these biomolecules resemble their native structures in vivo. It was recently revealed that the pharmacological activity of peptide drugs Enkephalins correlates with their gas-phase geometries, but not with the condensed-phase structures [4]. These peptide ligands adopt their biologically active geometry only when they dock into hydrophobic pocket of their receptors. Isolation of peptides in the gas phase mimics, to some extent, this hydrophobic environment as it removes strong intermolecular hydrophilic interactions.

In order to test this assumption, we applied conformer-selective multi-laser CIS to record conformer-selective vibrational spectra of individual conformers of three stereoisomers of the opioid pentapeptide [Ala₂-Leu₅]-Enk, isolated in the gas phase. Together with theoretical calculations, experimental spectra allowed for the identification of 3D structures of six conformers of three stereoisomers (3 for DD, 2 for DL, and 1 for LL). The obtained structures fully support the predictions made previously by FRET experiments in the gas phase. Furthermore, we found an evidence of the kinetically trapped, solution-like structure of one of the DL conformers.

Overall, these data confirm that for a large class of biomolecules there is high relevance of gas-phase studies for biology, as they reflect structural preferences of biologically active geometries.

References

1. Pereverzev A. Y., O. V. Boyarkin, Phys. Chem. Chem. Phys., 19, 3468 (2017).
2. Pereverzev A. Y., Cheng X., Nagornova N. S., Reese D. L., Steele R. P., Boyarkin O. V., J. Phys. Chem. A, 120, 5598 (2016).
3. Pereverzev A. Y., Kopysov V. N., Boyarkin O. V., Angew. Chem., 56, 15639 (2017).
4. Kopysov V. N., Boyarkin O. V., Angew. Chem., 55, 689 (2016).

MP-68 / THE USE OF PEPTIDE BONDS FOR COLD ION SPECTROSCOPY

Aleksandr Y. Pereverzev - Vladimir N. Kopysov - Oleg V. Boyarkin

EPFL, LCPM, Lausanne, Switzerland

Keywords: UV& IR spectroscopy, peptide bond, photodissociation

UV absorption by aromatic residues in cold, gas-phase peptides is extremely sensitive to their structure. Combined with vibrational spectroscopy (IR-UV double-resonance [1] and IR-IR-UV hole-burning [2, 3] techniques) it allows for accurate validation of calculated geometries, while together with high-resolution mass spectrometry (2D UV-MS [4]) it allows for an identification of isoforms of biomolecules. These photofragmentation techniques are applicable, however, only to those peptides that contain at least one aromatic residue.

Alternatively, peptide bonds can be used as chromophores. In solution, peptides and proteins exhibit strong $\pi \rightarrow \pi^*$ transition at ~ 190 nm and a weak $n \rightarrow \pi^*$ transition at ~ 215 nm. UV absorption of peptides in the gas phase, however, has never been characterized. We report the first gas-phase UV photofragmentation spectrum of a single peptide bond. The influence of vibrational excitation on the electronic spectrum of cold ions is explored and all-conformer IR-UV gain and conformer-selective IR-IR-UV hole-burning spectra of non-aromatic peptides with up to six amino acids are presented.

2D UV-MS spectroscopy of a model helical peptide reveals that the intensities of the $n \rightarrow \pi^*$ transitions of individual peptide bonds correlate with their local structural preferences. This type of transition is, usually, very weak because of small spatial overlap of the atomic non-bonding orbital of oxygen, n , and the molecular delocalized anti-bonding orbital of the peptide bond, π^* . We suggest that the intensity of the $n \rightarrow \pi^*$ transition depends on the degree of distortion of the O=C-N-H planar orientation due to the local, non-covalent interactions. Since 2D UV-MS spectroscopy is able to separate absorptions of individual peptide bonds, it is a promising structural probe for peptides and proteins.

References

1. Pereverzev A. Y., O. V. Boyarkin, Phys. Chem. Chem. Phys., 19, 3468 (2017).
2. Pereverzev A. Y., Cheng X., Nagornova N. S., Reese D. L., Steele R. P., Boyarkin O. V., J. Phys. Chem. A, 120, 5598 (2016).
3. Pereverzev A. Y., Kopysov V. N., Boyarkin O. V., Angew. Chem., 56, 15639 (2017).
4. Kopysov V. N., Makarov A., Boyarkin O. V., Anal. Chem., 87, 4607 (2015).

MP-69 / ACTION SPECTROSCOPY OF FLAVINS

Mark Stockett⁽¹⁾ - **Linda Giacomozzi**⁽¹⁾ - **Henning Zettergren**⁽¹⁾ - **Christina Kjær**⁽²⁾ - **Jeppe Langeland Knudsen**⁽²⁾ - **Lars Andersen**⁽²⁾ - **Steen Brøndsted Nielsen**⁽²⁾ - **James Bull**⁽³⁾ - **Eduardo Carrascosa**⁽³⁾ - **Evan Bieske**⁽³⁾

Stockholm University, Department of Physics, Stockholm, Sweden (1) - Aarhus University, Department of Physics and Astronomy, Aarhus, Denmark (2) - University of Melbourne, School of Chemistry, Melbourne, Australia (3)

Keywords: flavins, ion spectroscopy, dissociation, isomerization, luminescence

Introduction

Flavins are a ubiquitous family of redox cofactors which serve as blue light sensors in proteins and enzymes responsible for *e.g.* DNA repair, phototropism in plants and the perception of magnetic fields by some migratory birds. Little is known about how the optical properties of flavins are affected by micro-environmental perturbations, particularly in their reduced and/or deprotonated forms which act as resting or signaling states in some proteins.

Methods

Flavin ions are produced by electrospray. Photo-induced dissociation (PID) action spectroscopy was performed using a sector instrument and an electrostatic ion storage ring. Photo-induced isomerization action (PISA) spectroscopy was carried out using tandem ion mobility spectrometry (IMS). Laser induced fluorescence (LIF) was measured in a quadrupole ion trap (QIT).

Results

Using a multi-faceted approach employing state-of-the-art action spectroscopy experiments and concomitant calculations, we are investigating the photophysics and photochemistry of flavin ions in vacuo, with a focus on Flavin Adenine Dinucleotide (FAD). For FAD mono-anions, where the iso-alloxazine chromophore is in its fully oxidized, neutral form, the gas-phase absorption and fluorescence spectra are very similar to solution-phase measurements [1-3]. The gas-phase Stokes shift is considerably larger than previously reported for other molecular ions, implying significant vibronic effects [2]. Deprotonation of the chromophore leads to a significant red-shift relative to solution due to significant charge transfer character, which also gives a high susceptibility to micro-environmental perturbations [3]. Photo-induced intra-molecular proton transfer is found to be an important step in the photo-degradation of flavins [1,3]. Isomer-selective action spectroscopy has revealed deprotomer-dependent photochemistry [3].

Conclusions

Intrinsic spectroscopic and photochemical data are established. These gas-phase experiments form a baseline for quantitative comparison to theory and to more complex situations in solution or proteins. The photochemistry of flavins is driven by intra-molecular proton transfer. Time-resolved experiments are needed to disentangle excited-state and statistical processes.

Novel Aspect

First absorption and fluorescence spectra of flavin ions in vacuo. Photochemistry and micro-environmental sensitivity are strongly dependent on (de)protonation state.

References

1. Stockett M. H., Phys Chem Chem Phys, 19, 25829 (2017).
2. Giacomozzi L. et al., J Chem Phys, 148, 214309 (2018).
3. Bull J. N. et al, In preparation.

MPS-S04 - Hyphenated techniques

Chairs: Myeong Hee Moon, Michele Suman

MP-70 / SPEEDING UP GC-MS ANALYSIS – FROM FEW MINUTES DOWN TO REAL-TIME ANALYSIS

Alexander B. Fialkov - Uri Keshet - Tal Alon - Aviv Amirav

School of Chemistry, Tel Aviv University, Tel Aviv

Keywords: Real time analysis with separation

Introduction

Fast analysis requires speed optimization of the whole analysis including sample collection, introduction, separation, detection, identification and quantitation. We developed instruments and methods that each facilitates faster GC-MS analysis: Fast GC, Open Probe Fast GC and GC-MS with Cold EI.

Methods

The fast GC is based on a short capillary column inside thin walled metal tube that is resistively heated to provide sub one minute separation. Open Probe is based on a heated oven that is mounted on a low thermal mass fast GC which is open to room air with helium purge flow protection to prevent air penetration. Open Probe fast GC-MS operation is simple: touch the sample with a holder, insert it into the Open Probe and start running with 30 s separation and 50 s ready for next sample.

Preliminary data

We evaluated our low thermal mass Fast GC with Agilent 5975 GC-MS with standard EI and with GC-MS with Cold EI and demonstrated sub-one-minute full analysis cycle time. With the concept of low pressure chromatography we found that a mere factor of four loss of GC separation could bring a factor of 64 faster GC separations.

Open Probe Fast GC-MS provides sub-one-minute full analysis times, including sample preparation. Thus, Open Probe Fast GC-MS acts like the range of ambient desorption ionization techniques yet with separation, library identification at the isomer level, superior sensitivity and the use of low cost single quad MS of GC-MS.

Novel Aspect

New instruments and methods for low thermal mass fast GC and Open Probe Fast GC-MS for sub one minute fast GC-MS and real time analysis with separation

MP-71 / DEVELOPMENT OF COMPREHENSIVE STEROID ANALYSIS METHODS USING GCxGC-TOFMS

Viatcheslav Artaev⁽¹⁾ - ***David Alonso***⁽²⁾ - ***Jonelle Shiel***⁽¹⁾ - ***Georgy Tikhonov***⁽¹⁾ - ***Albert Lebedev***⁽³⁾

LECO Corp, R&D, St. Joseph (1) - LECO Corporation, Application Dept., St. Joseph (2) - Moscow State University, Chemistry Dept., Moscow (3)

Keywords: TOFMS, HRMS, GCxGC, Steroids, Clinical

Introduction

Identification and quantitation steroid metabolites in biological samples are essential for screening for various hormonal disorders. Many of these metabolites are closely related isomers and cannot be easily separated by LC or one dimensional GC due to their chemical and structural similarity. Comprehensive GCxGC coupled to a high resolution TOFMS provides enhanced chromatographic separation, better detection and identification of the analytes of interest.

Methods

All data was obtained using the high resolution multiple reflecting geometry time-of-flight (TOF) mass spectrometer equipped with GCxGC with quad-jet, liquid nitrogen cooled thermal modulator and secondary column oven. We have tested various GC columns sets and oven temperature programs during GC method development. Steroids standards were obtained from Steraloids (USA) and Fisher Scientific (USA). Samples were derivatized using MSTFA method.

Results

The steroid standards were analyzed individually and from that runs deconvoluted and curated mass spectra were obtained and included in the custom accurate mass library to assist in analytes identification in the complex mixtures. We present a workflow and analysis method that enables detection of 33 steroids from different classes in urine. Enhanced peak capacity provided by GCxGC allows separation of all analytes of interest, including structural isomers, which cannot be resolved or identified by other techniques. The high mass accuracy data along with automatic deconvolution algorithms provide reliable mass spectral information allowing accurate peak assignment. The method was transferred on to another analytical setup – the low resolution (LR) TOFMS coupled to GCxGC. Use of GCxGC-LRTOFMS represents a more practical approach of implementing the method for steroid analysis, while the HRTOFMS results can be applied for reliable analyte assignment by using HRTOFMS RI data and mass spectral match.

Conclusions

The method is developed for comprehensive analysis of steroids in urine by GCxGC-HRTOFMS as well as GCxGC-LRTOFMS. The method allows taking full advantage of increased separation capacity of GCxGC and mass accuracy and sensitivity of the TOFMS. The reproducible sample preparation (derivatization, cleanup, etc) was the main challenge, and the future work will be focused on automation of the sample preparation process.

Novel Aspect

33 steroids were simultaneously analyzed in complex biological matrix using GCxGC separation coupled with low and high resolution TOFMS with automatic data processing.

MP-72 / USING ELECTROKINETIC INJECTION TO INCREASE THROUGHPUT AND IMPROVE SENSITIVITY IN THE DETECTION OF PROTEINS BY CE-MS

Jim Thorn⁽¹⁾ - *Stephen Lock*⁽¹⁾ - *Walter Belloni*⁽²⁾ - *Christopher Löβner*⁽²⁾

SCIEX, Warrington, United Kingdom (1) - SCIEX, Darmstadt, Germany (2)

Keywords: CE-MS, CESI, electrokinetic, sensitivity

Introduction

Capillary electrophoresis (CE) is an orthogonal technique to LC separating compounds based on their charge. The properties of CE enable the reduction and often elimination of carryover and wall absorption which effects peak resolution and sensitivity of LC. In this work we will describe how a CESI-MS method has been developed to detect intact proteins between 15 -25,000 amu. CESI [the integration of capillary electrophoresis (CE) and electrospray ionization (ESI) into a single process in a single device] is now enabling the easy connection of CE to mass spectrometers.

Methods

Protein standards were prepared in a variety of different solvents and injected by either pressure or electro-kinetically onto a neutrally coated capillary under various conditions. The CE separation used MS amenable background electrolyte (BGE) consisting of a mixture of Acetic acid, Acetonitrile and Water. Proteins were detected by a MS system in MRM mode (using a capillary 30 μm ID, 91 cm long) which was fitted with a NanoSpray® III source which was used in either full scan or single ion monitoring mode with an ionspray voltage of 1600 - 1800 V which was optimized to the BGE being used.

Results

The CESI-MS method was developed with the capability of detecting intact protein standards at <10 ng/mL. This study showed that electro kinetic injection was >40 fold more sensitive than isotachopheresis (ITP)

injection technique and is also more sensitive than a traditional LCMS approach and followed a similar trend as shown previously for neuropeptides. The response obtained was linear over the 2-3 orders tests with very low levels of carryover observed. Sensitivity depended on how well the protein ionized and one of the major factors shown to affect sensitivity on the electro kinetic injection was the levels of acetonitrile in the sample.

Novel Aspect

Increasing sensitivity for bioanalysis by using electrokinetic injection techniques

MP-73 / MULTIVARIATE STUDY OF MATRIX EFFECT IN THE ANALYSIS OF PHYTOESTROGENS IN SOY-FOOD BY LC-ESI-MS/MS

Barbara Benedetti - Marina Di Carro - Emanuele Magi

University of Genova, Department of chemistry and industrial chemistry, Genova

Keywords: liquid chromatography-tandem mass spectrometry, matrix effect, experimental design, QuEChERS, phytoestrogens

Introduction

Matrix effect is an important aspect to consider when analyzing extracts from complex food matrices by LC-ESI-MS [1]. In fact, if no labeled standards are available, accurate quantitation is a difficult task. The work deals with the evaluation of matrix effect involved in the analysis of phytoestrogens in soy-based food. QuEChERS [2] were used as extraction/purification technique and the multivariate approach of experimental design was used to optimize the procedure.

Methods

The phytoestrogens daidzein (DAID), genistein (GEN), coumestrol (COUM), formononetin (FORM) and biochanin A (BIOCH) were extracted from soy-based burgers using the QuEChERS methodology: the steps involved are extraction by acetonitrile-water, phase separation by salts addition and clean-up. Two subsequent experimental designs were employed to study the factors influencing the procedure. LC-MS/MS in multiple reaction monitoring mode was used for the analysis.

Results

Matrix effect (ME) was assessed by comparing the analytes' signals in neat standards and in soy-burger extracts. A screening design (Plackett-Burman)[3] was employed to select the important variables involved in the QuEChERS procedure and affecting ME (response). No factors were significant for DAID and GEN; in fact, they were determined in diluted extracts, where negligible ME was observed. As for COUM, FORM and BIOCH, the factors considered statistically significant ($p < 0.05$) were: sample amount, PSA (Primary Secondary Amine) sorbent and FLORISIL sorbent (used in the clean-up stage). A subsequent Box-Behnken design [4] allowed to obtain the response surfaces (models of the response as a function of the variables). The models for COUM and FORM were subjected to ANOVA and showed good R^2 coefficients (0.85 and 0.78 respectively). Conversely, the model for BIOCH was not significant. The optimal conditions to minimize ME (ion suppression $< 20\%$) were: 200 mg of sample (for 10 mL of solvent) and clean-up of 2 mL of extract with 40-50 mg of PSA and 80-90 mg of FLORISIL.

Conclusions

Thanks to the approach of experimental design and response surface methodology, it was possible to study the problem of matrix effect in LC-ESI-MS analysis of complex matrices (soy-based burgers). The QuEChERS procedure, including a clean-up step, was optimized with a limited number of experiments, despite the

several factors involved. The method is cheap, rapid, does not require specific equipment and provides a valid alternative to the use of labeled internal standards.

Novel Aspect

Application of the multivariate approach to optimize matrix effect in the analysis of a complex and novel matrix (soy-based burger) without the use of internal standards.

References

1. B. K. Matuszewski, M. L. Constanzer, C. M. Chavez-Eng, *Anal. Chem.* 75 (2003) 3019–3030.
2. M. Anastassiades, S. J. Lehotay, D. Stajnbaher, F. J. Schenck, *J. AOAC Int.* 86 (2003) 412–31.
3. R. L. Plackett, J. P. Burman, *Biometrika Trust.* 33 (1946) 305–325.
4. G. E. P. Box, D. W. Behnken, *Technometrics.* 2 (1960) 455–475.

MP-74 / APPLICATION OF SPME GC-MS TO EXPLORE POSSIBLE DIFFERENCES IN FREE-LACTOSE MILKS FROM THE MARKET

Riccardo Bottioli - Eugenio Aprea - Emanuela Betta - Flavia Gasperi

Edmund Mach Foundation, Food Quality and Nutrition Department, San Michele all'Adige

Keywords: lactose-hydrolyzed milk, SPME GC-MS, volatile organic compounds

Introduction

Positive effects of milk on the human health are well-known. However, 70% of the global population is lactose intolerant [1]. For this reason lactose-free products, obtained adding lactase, have been released on the market [2]. According to the step in which lactase is added, downside effects can arise over time, due to its proteolytic side activity [3]. The present study aim to explore the potential of SMPE GC-MS to catch VOCs profiling of commercial lactose-free milk.

Methods

UHT lactose-free milks from the local market, three freshly produced (FM1, FM2, FM3) and two after 180 days of storage (SM1 and SM2) were analyzed by SPME GC-MS according to Bergamaschi et al. (2015) [4]. Compounds were identified using the NIST-98/Wiley library and calculated retention indices (RI). Samples were compared by one-way ANOVA and Tukey post hoc test, when necessary, using the STATISTICA software (Dell Software Inc., Palo Alto, CA).

Results

Fifty-nine volatiles were identified and 64% of them were present in all the samples, suggesting differences in the VOCs profile. Many of those were already reported in other studies focused on lactose-free milk [3, 5, 6]. Fresh samples (FM1, FM2 and FM3) differed significantly for 7 compounds and 6 of them were higher in FM3. At the moment, we cannot speculate on the observed differences but the results suggested that the SPME GC-MS can catch the differences among milk batches. Fresh and stored milk differed significantly for 6 compounds, 5 of them having higher intensity after 180 days of storage (SM1 and SM2). Several of the VOCs that differed among the samples were methyl-ketones, well-known for their contribution to the off-flavor in UHT milk [7]. Some, such as 2-undecanone and 2-tridecanone, have long carbon chain and can be responsible for the heated flavor of milk [7]. The formation of methyl-ketones can also be addressed to the Maillard reaction, whose proceeding can be referred to the higher levels of these compounds found in SM1 and SM2.

Conclusions

The study indicated that the SPME GC-MS approach was able to catch differences in volatiles profiling of the milk tested. Conservation had a relevant impact and most of the changes were related to compound

intensity. As the proteolytic side activities of lactase may play a role, it will be verified further. A shelf-life study is envisaged to verify the evolution of the VOCs profile in the samples and its contribution to the products sensory profile.

Novel Aspect

This study is an attempt to verify the suitability of the SPME GC-MS for tracking the evolution of the off-flavours formation in lactose-hydrolysed milk during storage.

References

1. Swallow D. M., Annual review of genetics, 37, 197-219 (2003).
2. Jelen P., Tossavainen O., Australian Journal of Dairy Technology, 58(2), 161(2003).
3. Troise A. D., Bandini E., De Donno R., Meijer G., Trezzi M., Fogliano, V., Food Research International, 89, 514-525 (2016).
4. Bergamaschi M., Aprea E., Betta E., Biasioli F., Cipolat-Gotet C., Cecchinato A., Bittante G., Gasperi F., Journal of Dairy science, 98(4), 2183-2196 (2015).
5. Jansson T., Jensen S., Eggers N., Clausen M. R., Larsen L. B., Ray C., Sundgren A., Andersen H. J., Bertram H. C., Dairy Science & Technology, 94(4), 311-325 (2014).
6. Jensen S., Jansson T., Eggers N., Clausen M. R., Larsen L. B., Jensen H. B., Ray C., Sundgren A., Andersen H. J., Bertram H. C., European Food Research and Technology, 240(6), 1247-1257 (2015).
7. Contarini G., Povolito M., Leardi R., Toppino, P. M., Journal of Agricultural and Food Chemistry, 45(8), 3171-3177(1997).

MP-75 / POST-COLUMN IN-SOURCE DERIVATIZATION IN LCMS: A TOOL FOR NATURAL PRODUCTS CHARACTERIZATION AND METABOLOMICS

Adriana Cirigliano⁽¹⁾ - **Gabriela M Cabrera**⁽²⁾

DQO, FCEN, Universidad de Buenos Aires, DQO, FCEN, Universidad de Buenos Aires, Buenos Aires (1) - DQO, FCEN, Universidad de Buenos Aires, UMYMFOR-DQO,UBA-CONICET, Buenos Aires (2)

Keywords: Post-column derivatization, LC-MS, in-source derivatization

Introduction

Metabolite characterization is crucial for many aspects of basic research, such as de-replication of crude extracts in natural products chemistry or tentative identification in untargeted metabolomics. This information makes the subsequent decision-making process easier and more efficient. Derivatization with hydrazines is a known method in this respect, which allows the characterization of carbonyl functional groups [1-2].

Methods

Liquid Chromatography-Mass Spectrometry was performed using Electrospray or Atmospheric Pressure Chemical Ionization as ionization techniques. The post-column addition of the reactants (hydrazines) was done by means of a syringe pump via a T-junction before entrance into the ion source. Experimental conditions were optimized. LCMS runs without addition were also performed for comparison.

Results

Different compounds with at least one carbonyl group in its structure were analysed by this method. Triterpenoids like cycloartanes, steroids, cardenolides and other terpenoids were included in this study. All the analysed samples showed a similar behavior. The spectra of compounds with aldehyde groups exhibited as main signals those corresponding to the hydrazones, which are the expected products for the reaction between aldehydes and hydrazines. On the contrary, compounds with conjugated ketones did not react with hydrazines and did not show additional signals, and compounds having non-conjugated ketones

presented in their spectra mainly those signals corresponding to the product of the nucleophilic addition adduct of hydrazine to the carbonyl group.

Conclusions

The in situ derivatization of carbonyl compounds in the ionization source was achieved, and applied in high-performance liquid chromatographic analysis. The typical derivatives of carbonyl compounds with hydrazines were obtained. Other reactions are being evaluated to further extend this method to other classes of derivatives.

Novel Aspect

The use of post-column in-source derivatization in LCMS allows the characterization of compounds with carbonyl groups in their structures in a quick and easy way.

References

1. Lattová E., Perreault H., Mass Spectrometry Reviews, 32, 366-385 (2013).
2. Yasuhara A., Tanaka Y., Makishima M., Suzuki S., Shibamoto T., Journal of Chromatography & Separation Techniques, 2, 1-5 (2011).

MP-76 / CONJUGATION AMONG NUCLEOTIDES AND CYANOCOBALAMIN FUNCTIONALIZED WITH CISPLATIN: A STUDY BY LIQUID CHROMATOGRAPHY COUPLED WITH ESI AND MULTISTAGE MS

Giovanni Ventura - Cosima Damiana Calvano - Andrea Viola - Valeria Cinquepalmi - Francesco Palmisano - Tommaso Cataldi

University of Bari, chemistry department, Bari, Italy

Keywords: LC-MS, cisplatin, cyanocobalamin, RNA, prodrugs

Introduction

Cisplatin has been largely used in antineoplastic therapy [1]; diverse drug delivery systems have been studied in order to reduce its serious side effects [2]. DNA is general accepted as the critical target for platinum-based compounds, nonetheless adducts between cisplatin and proteins [3], lipids [4] and RNA [5] might contribute to their effectiveness. In this study, interaction between cyanocobalamin-cisplatin (CNCbl-CP) prodrug and RNA bases has been studied via LC-MS.

Methods

After the reaction between CNCbl and CP, a complex with empirical formula $[C_{63}H_{94}ClCoN_{16}O_{14}Ppt]^+$ is generated [6]; its reactivity with in vitro incubation of individual adenosine monophosphate (AMP), guanosine monophosphate (GMP), cytosine monophosphate (CMP) and uridine monophosphate (UMP) in aqueous solution at 37°C was studied by comparing results obtained using two different HPLC columns coupled to electrospray ionization and a linear ion trap (ESI-LIT) mass spectrometer.

Results

The use of innovative drug delivery systems able to increase selectivity against cancer cells is a common strategy to reduce cisplatin toxicity; CNCbl functionalized with cisplatin is an example of prodrug. The formation of adducts between CNCbl-CP and RNA purine and/or pyrimidine bases was studied following two approaches: reaction was performed first on individual NMPs and later with competitive reactions in simultaneous presence of four NMPs. An appropriate gradient elution program based on methanol was developed to attain a separation both on a conventional and an amide-embedded ODS stationary phase; those two columns may provide complementary selectivity towards platinum adducts. All the NMPs react

with CNCbl-CP by substitution of the last chlorine ligand on the cisplatin. The comparison between obtained and simulated isotopic patterns together with MS/MS analyses confirmed the presence of platinum in the reaction products.

Conclusions

From the obtained data it was clear that AMP and GMP are the most reactive nucleotides due to the presence of purine rings; three isobaric chemical species have been identified at m/z 965.8 for AMP, while a single protonated adduct at m/z 973.8 was generated for GMP. Also, adducts with CMP and UMP have been identified at m/z 953.8 and m/z 954.3, respectively. Plausible structures have been proposed by using tandem (MS/MS) and MS3 spectra.

Novel Aspect

The interaction between CNCbl-CP and RNA bases has been proved by using HPLC-MS; the examination of fragmentation patterns was found very useful to gain insight into detailed molecular information.

References

1. D. Wang, S. J. Lippard, *Nat. Rev. Drug Discov.*, 4, 307–320, (2005).
2. T. C. Johnstone, K. Suntharalingam, S. J. Lippard, *Chem. Rev.*, 116, 3436–3486, (2016).
3. C. G. Hartinger, Y. O. Tsybin, J. Fuchser, P. J. Dyson, *Inorg. Chem.*, 47, 17–19, (2008).
4. M. Jensen, W. Nerdal, *Eur. J. Pharm. Sci.*, 34, 140–148, (2008).
5. A. A. Hostetter, M. F. Osborn, V. J. Deroose, *ACS Chem. Biol.*, 7, 218–225, (2012).
6. G. Ventura, F. Arnesano, C. D. Calvano, F. Palmisano, T. R. I. Cataldi, *RSC Adv.*, 7, 53658–53666, (2017).

MP-77 / ANALYTICAL METHOD FOR QUANTIFICATION OF BIXIN IN EXTRACT OF ANNATTO USING V-EASI(-)-LC-MS A METHOD SOFTER THAN ESI-LC-MS

Luan Felipe Campos Oliveira - David Ulisses Tega - Mirela B. Coelho - Damila Rodrigues de Moraes - Mariana Baptista - Marcos N. Eberlin

UNICAMP, Chemistry Institute, Campinas, Brazil

Keywords: Analytical methods, mass spectrometry, V-EASI, LC-MS, Bixin

Introduction

Ventury Easy Ambient Sonic-Spray ionization (V-EASI) is a ionization source softer than electrospray (ESI) derived from Sonic Spray (SSI) [1]. Until now it has been little explored, used only for qualitative analyses [2-4]. This study is a proof of concept of an LC-MS system using the V-EASI source applied to develop a quantitative method of bixin, a natural dye from annatto sensitive to high voltage [5].

Methods

The quantification was performed in LC-MS system using a Waters 510 chromatographic pump coupled with V-EASI source, and the quadrupole-mass spectrometer (LCMS Shimadzu 2010), where was used the column C18 in the measurements. Methanol, bixin standard, lauric acid (internal standard), acetic acid (modifier) were used. For the validation, the guidance of validation from FDA was employed [6].

Results

The method was developed and optimized using the SIM mode for the ion of m/z 199 (lauric acid) and ion of m/z 393 (bixin) by V-EASI(-). We have found with the optimization the best analysis condition, as follow: 6 Bar of source, MS analyzer and block at the temperature of 250°C, detector voltage of 1.50 kV, chromatographic run with 100% of methanol (0.05% acetic acid, v/v) where internal standard eluting at 3.8 min and bixin in 4.4 min. The validation data showed a linear range of 0.5 to 10 ppm, limit of detection and

quantification of 0.02 and 0.07 ppm, respectively. The imprecision interday and intraday were 4.6% and 8.4%, respectively. The method also was considered accurate in an agreement with the official method by HPLC-DAD [7], it provides statistically equal results and has been shown to be a selective method with no matrix effect.

Conclusions

We were able to do the validation of a new V-EASI(-)-LC-MS method for bixin quantification, which presented similar results in comparison with the official method [6]. The method opens the perspective of a new system aiming the quantification of sensitive compounds and may be a new possibility for omics sciences.

Novel Aspect

A powerful coupling among mass spectrometer, liquid chromatography and V-EASI source for the bixin quantification.

References

1. Santos, V. G., T. Regiani, F. F. Dias, W. Romao, J. L. Jara, C. F. Klitzke, F. Coelho & M. N. Eberlin (2011) Venturi easy ambient sonic-spray ionization. *Anal Chem*, 83, 1375-80.
2. Haddad, R., T. Regiani, C. F. Klitzke, G. B. Sanvido, Y. E. Corilo, D. V. Augusti, V. M. D. Pasa, R. C. C. Pereira, W. Romão, B. G. Vaz, R. Augusti & M. N. Eberlin (2012) Gasoline, Kerosene, and Diesel Fingerprinting via Polar Markers. *Energy & Fuels*, 26, 3542-3547.
3. Porcari, A. M., G. D. Fernandes, D. Barrera-Arellano, M. N. Eberlin & R. M. Alberici (2016) Food quality and authenticity screening via easy ambient sonic-spray ionization mass spectrometry. *Analyst*, 141, 1172-84.
4. Schafer, K. C., J. Balog, T. Szaniszló, D. Szalay, G. Mezey, J. Denes, L. Bogner, M. Oertel & Z. Takats. 2011. Real time analysis of brain tissue by direct combination of ultrasonic surgical aspiration and sonic spray mass spectrometry. In *Anal Chem*, 7729-35.
5. Scotter, M. (2009) The chemistry and analysis of annatto food colouring: a review. *Food Additives & Contaminants: Part A*, 26, 1123-1145.
6. FDA Guidance for Industry, Q2B Validation of Analytical Procedures Methodology, 1996.
7. Combined Compendium of Food Additive Specifications Joint FAO/WHO Expert Committee on Food Additives 2006.

MP-78 / THIN LAYER CHROMATOGRAPHY COMBINED WITH FLAME-INDUCED ATMOSPHERIC PRESSURE CHEMICAL IONIZATION MASS SPECTROMETRY (FAPCI/MS) FOR VOLATILE AND SEMI-VOLATILE COMPOUND ANALYSIS

Sychyi Cheng - Jentaie Shiea

National Sun Yat-sen University, Department of Chemistry, Kaohsiung

Keywords: thin layer chromatography, desorption flame-induced atmospheric pressure chemical ionization

Introduction

The combustion of hydrocarbons in a flame produces reactive ion species. In addition, the high temperature generated in a flame allows thermal desorption of small compounds on surfaces. A flame is therefore directed toward the sample to characterize analytes on surfaces, so called desorption flame-induced atmospheric pressure chemical ionization (DFAPCI). In this study, thin layer chromatography (TLC) was combined with DFAPCI/MS to analyze a mixture.

Methods

The developed TLC plates were set on a XY stage and moved by a robotic platform for DFAPCI/MS analysis.

A combustion head fabricated by arranging two concentric stainless steel tubes was utilized to generate a micro-flame. The combustion head was set at an angle of 45° with respect to the TLC plate surface. Volatile and semi-volatile compounds such as caffeine, drug molecules, and plant extracts were analyzed by TLC-DFAPCI/MS.

Results

A light-brown line with 1-1.5 mm width was produced on the gel bed after analysis. This line clearly revealed the surface bearing a steam of heated jet from an oxyacetylene flame. The surface temperature of a TLC plate was measured to be 300 ± 10 °C as the plate passing through the desorption/ionization (DI) region. A mixture containing nicotinamide, 2-phenylacetamide, and dibenzylamine was separated by a silica gel TLC plate for the subsequent DFAPCI/MS detection. Positive ions of nicotinamide (m/z 123, $[M+H]^+$), 2-phenylacetamide (m/z 136, $[M+H]^+$), and dibenzylamine (m/z 196, $[M-H]^+$) were detected. The limit of detections (LODs) of TLC-DFAPCI/MS for the analysis of three standards were tested and concluded to be 50 ng/spot for nicotinamide, 25 ng/spot for 2-phenylacetamide, and 5 ng/spot for dibenzylamine. Samples such as amine and amide standards, drug molecules, and aromatherapy oil were successfully characterized by TLC-DFAPCI/MS.

Conclusions

Compared with other TLC-AMS approaches, TLC-DFAPCI/MS shows advantages including cheap, simple instrumental setup, and easy to be fabricated. Since the analytes were vaporized by a heat, TLC-DFAPCI/MS was suitable to characterize volatile, semi-volatile, and thermally stable compounds. We are combining laser desorption (LD) with FAPCI/MS, which will be capable of detecting non-volatile substances.

Novel Aspect

Thin layer chromatography (TLC) has combined with desorption flame-induced atmospheric pressure chemical ionization (DFAPCI) to analyze a mixture.

References

1. Cheng S. C., Chen S. H., Shiea J., Mass Spectrometry (Tokyo), 6, S0065 (2017).
2. Cheng S. C., Wang C. H., Shiea J., Anal. Chem., 88, 5159 (2016).
3. Cheng S. C., Chen Y. T., Jhang S. S., Shiea J., Rapid Commun. Mass Spectrom., 30, 890 (2016).

MP-79 / EVALUATION OF GAS CHROMATOGRAPHY COUPLED TO ELECTROSPRAY-ATMOSPHERIC PRESSURE PHOTON IONIZATION ORBITRAP MASS SPECTROMETRY AS AN ADVANCED STRATEGY FOR POLYAROMATIC HYDROCARBON ANALYSIS

Yunju Cho⁽¹⁾ - Eunji Cho⁽²⁾ - Sunghwan Kim⁽¹⁾

Green-Nano Material Research Center, Kyungpook National University, Chemistry, Daegu (1) - Kyungpook National University, Chemistry, Daegu (2)

Keywords: Orbitrap, Gas chromatography, Mass spectrometry, High resolution, polyaromatic hydrocarbons

Introduction

Heteroatom containing PAHs (HPAHs) play an important role in petroleomics due to their toxicity, stability in links of the catalyst development during oil refining process. However, general GC-EI/MS is not suitable for polar PAHs due to ionization process. [1] Also, no single ionization technique has capability to ionize different variety of compounds. [2] In this study, GCESI-APPI MS provides a technique with high sensitivity and accurate m/z values of the HPAHs.

Methods

About thirty PAHs including nitrogen, sulfur and oxygen containing compounds were used for the

establishment of GC ESI-APPI orbitrap MS conditions. GC was performed using an Agilent 7890A equipped with in-house-heated column transfer line that was heated by a programmable temperature controller with a temperature sensor. [3, 4] Mass spectrometer was orbitrap Q-Exactive (ThermoFisher Scientific Inc.) equipped with an ESI source and APPI lamp.

Results

Depending on the analytes, GC has been coupled to various detectors, *e.g.* MS, FID and NPD. GC-MS is the most commonly employed in EI. But EI may be lost the parent ion of the analytes because of extensive fragmentation. FID is the presence of heteroatoms in molecular decreases the detector's response. To solve this problem, element-selective detectors, like NPD, have been used. But NPD has limited signal stability due to mechanism of detection that is based on the bead. [1] To overcome such concerns, orbitrap MS with high sensitivity and resolution were applied to the HPAHs combining ESI which is soft ionization causing minimal fragmentations and APPI that it has been used to study PAHs in oils [5] hyphenated GC. The several ESI solvents were compared [6] and acetone showed the best ionization efficiency among the tested solvents. Also it is observed that the results using GC ESI-APPI MS for HPAHs is higher sensitivity than using GC-MS with single ESI source. Therefore, the GCESI-APPI MS is suited for simultaneous analysis of HPAHs as an advanced ionization.

Conclusions

The aim of this research was to evaluate an instrument to couple GC with high-resolution orbitrap MS using an ESI source and APPI lamp. The studies demonstrated that the GC ESI-APPI MS has more elements-selectivity, sensitivity and accuracy even better than obtained with traditional GCEI/MS for HPAHs. Therefore, we expect that this method should be useful in the selective analysis of complex mixtures as nitrogen-rich oil, spilled oil, metabolite and humic substances.

Novel Aspect

GC coupled to high-resolution orbitrap MS with ESI-APPI was used to evaluate heteroatom containing PAHs analysis at the molecular level compared with GC ESI MS and GC EI MS.

References

1. Grzegorz B., Patrycja M., Andre F., Andrzej P., *Journal of Separation Science*, 40, 1301 (2017).
2. Yunju C., Jang M. J., Matthias W., Justin E. B., Jeong G. N., Nam S. R., Sunghwan K., *Energy & Fuels*, 27, 1830 (2013).
3. Eunju C., Eun S. J., Sangwon C., Jaeick L., *Analytica Chimica Acta*, 964, 123 (2017).
4. Eun S. J., Eunju C., Sangwon C., Sunghwan K., Han B. O., Oh S. K., Jaeick L., *Analytical Chemistry*, 89, 12284 (2017).
5. Yunju C., Justin E. B., Manhoi H., Joonhee L. Byungjoo K., Sunghwan K., *Energy & Fuels*, 31, 7874 (2017).
6. Mirabelli M. F., Zenobi R., *Analytical Chemistry*, 90, 5015(2018).

MP-80 / ONLINE COUPLING OF NON-DENATURING CHROMATOGRAPHY TO NATIVE MASS SPECTROMETRY IS NO MORE AN IMPOSSIBLE MARRIAGE

Anthony Ehkirch - Oscar Hernandez-Alba - Thomas Botzanowski - Sarah Cianf erani

LSMBO/IPHC, University of Strasbourg/CNRS UMR7178, Strasbourg

Keywords: native mass spectrometry, ion mobility, LC-MS, therapeutic protein

Introduction

Mass spectrometry performed in non-denaturing conditions (native MS or IM-MS with ion mobility) has gained interest for the qualitative and quantitative characterization of intact proteins. However, online hyphenation of non-denaturing chromatographic (ndLC) techniques to native MS/IM-MS is still not

straightforward, which hampers its routine use in high throughput automated environments.

Methods

The 1DLC-MS system consists of an Acquity H-Class coupled to a Synapt G2 HDMS (Waters). Several Acquity BEH SEC (Waters) columns were used for 1D SEC-native MS. For 2D LC-MS, an AcquityI-Class was added to the 1D setup. For the 1st dimension, the column was a MAbPac HIC-10 (Thermo Scientific) or a Yarra SEC-X300 (Phenomenex); the 2nd SEC dimension was an AdvanceBio SEC (Agilent Technologies).

Results

We present here 1D and 2D LC-MS setups for online coupling of ndLC to native MS for the intact protein analysis. Native MS is mostly performed in ammonium acetate (AcONH₄) [1] while ndLC techniques (*e.g.* size exclusion SEC or hydrophobic interaction HIC chromatography) require large amounts of non-volatile salts [2]. We present here the benefits of 1D SEC-native MS in AcONH₄ for automated online buffer exchange and protein aggregation studies [3]. As AcONH₄ is not optimal for HIC/SEC, a multidimensional approach combining HIC/SEC in the 1st dimension and SEC in the 2nd one (HICxSEC or SECxSEC) was developed. The 1st ndLC dimension (HIC/SEC) aims at achieving optimal protein separation, while the 2nd ndLC (SEC) serves for fast desalting. Benefits of HIC/SECxSEC-native MS/IM-MS are shown for different types of proteins, including monoclonal antibody-based (mAb) products and amyloid proteins.

Conclusions

Our results demonstrate the possibility to online couple ndLC to native MS/IM-MS. Benefits of ndLCxLC-native MS/IM-MS setups provided comprehensive and streamlined characterization of different types of proteins [4]. The synergic online coupling HIC/SEC to native MS/IM-MS is envisioned to definitely push native MS approaches at the forefront of intact protein characterization.

Novel Aspect

Online hyphenation of non-denaturing LC to native IM-MS; 1D SEC-native MS/IM-MS; 2D HIC/SECxSEC-native MS/IM-MS

References

1. vanDuijn E., J Am Soc Mass Spectrom., 21, 971 (2010).
2. Goyon A, D'Atri V, Colas O, J Chromatogr B, 1065-1066, 35 (2017).
3. Ehkirch A., Hernandez-Alba O., Colas O., J Chromatogr B, under press (2018).
4. Ehkirch A, D'Atri V, Rouviere F, Anal Chem., 90, 1578 (2018).

MP-81 / IDENTIFICATION AND TIME VARIATION ANALYSIS OF PERFUME BY USING GCxGC/HRTOFMS WITH EI/FI/PI

Robert B. Cody⁽¹⁾ - **Michael C. Long**⁽²⁾ - **A. John Dane**⁽¹⁾ - **Koji Okuda**⁽¹⁾

JEOL USA, Inc., Mass Spectrometry, Peabody, MA (1) - Biopalette, Inc., Independent, Troy, NY (2)

Keywords: Perfume analysis, Comprehensive GC, High resolution, High mass accuracy, Soft ionization

Introduction

GCxGC-TOFMS is a powerful tool for measuring complex samples. Using EI with this method allows for NIST searches, but these searches may result in mis-assignments without additional exact mass information for the molecular ion as generated through soft ionization techniques like FI or PI. In this study, we measured perfume oil by using GCxGC-HRTOFMS with EI, FI and PI.

Methods

All samples were measured by using a high resolution TOFMS and GCxGC system. A pure perfume oil was qualitatively measured using EI, FI and PI. For the time studies, 0.1 μL of pure oil for the same perfume was placed in a 22mL vial and heated at 35°C for 0 hour, 0.5 hour, 1 hour and 3 hours. Afterwards, the headspace was sampled by using a SPME fiber and then measured with GCxGC FI.

Results

The 2D chromatograms for EI, FI and PI showed the same retention time (R. T.) patterns. Additionally, the FI and PI data showed stronger molecular ion signals than the EI data. As an example, the peak detected at R. T. 25.70 minutes(1st column) and 1.55 seconds(2nd column) showed relatively high fragment ion intensities for the EI measurements. The NIST search for this mass spectrum showed several compounds with a match factor over 900 but with different chemical compositions. In order to narrow down the possibilities, it was necessary to measure the molecular ion using FI and PI. The results for these techniques showed an accurate m/z value within 1 mDa of the calculated mass for $\text{C}_{12}\text{H}_{20}\text{O}_2$. This estimated composition narrowed the EI library search candidate to linalyl acetate.

Conclusions

The combination of EI fragmentation information with the accurate mass molecular ion information from FI and PI provided a means for narrowing down the NIST library search to help identify the compounds measured. We will present the qualitative results for the other compounds in the sample along with the perfume time variation studies.

Novel Aspect

Qualitative and quantitative analysis of perfume by GCxGC–HRTOFMS with EI, FI and PI ionization methods.

MP-82 / SECONDS-PER-SAMPLE: NOVEL ACOUSTIC HYPHENATED APPROACHES TO ENABLE TRULY LABEL-FREE ULTRA HIGH-THROUGHPUT AMBIENT MASS SPECTROMETRY ADVANCING DRUG DISCOVERY

Sammy Datwani - Lucien Ghislain

Labcyte, Inc - Chemistry, San Jose

Keywords: Acoustic droplet ejection (ADE), Acoustic mist ionization (AMI), Open port probe (OPP)

Introduction

A truly label-free ultra-high throughput screening technology is the “holy grail” platform for drug discovery. We have integrated acoustic droplet ejection (ADE) technology that couples ambient mass spectrometry with two novel approaches enabling low sample consumption (nanoliters), low cost and ultra-high-throughput.

Methods

We first highlight the principles of ADE and key technologies that enable robust acoustic liquid handling operations. The first modality utilizes direct electro spray ionization (ESI) whereby ADE delivers up to 1,500 spray events per second from an assay well in a microplate into the MS via a heated transfer optic [1]. The second modality explores an exciting method of sample injection to the MS by integrating ADE with 2.5 nL droplets into an open-port probe (OPP) sampling interface [2].

Results

The first modality, acoustic mist ionization (AMI), integrates ADE with a Waters Xevo G2-XS QToF mass spectrometer. The utility of this integrated system as an ultra-high throughput platform in this format was demonstrated on a biochemical screen of more than 250,000 compounds identifying new and novel

inhibitors of a relevant drug target. The sample assay was able to tolerate the changes to buffer conditions without impacting enzyme activity or kinetics. A data processing algorithm was designed to deal with large sample batches.

In the second modality the ADE-OPP-ESI-MS system operates with a continuous-flow of carrier solvent (MeOH) for direct sample dilution that also actively cleans the entire flow system, mitigating carry over. "Classic" ESI ion source used in this setup enables detection of a broad range of analytes including small drug molecules, antibody standards, peptides (MW>1500) and proteins. High sensitivity (attomol loading for small molecules, or sub-femtomol for intact antibody) and reproducibility (<8% CV) was demonstrated without internal standard normalization.

Conclusions

We highlight ultra-high throughput ambient mass spectrometry hyphenated approaches with an acoustic sample introduction interface - a MS based plate reader that enables broad applications in drug discovery. Next steps on the commercialization perspectives are addressed with these two modalities and the need to integrate the MS spectra "big data" collection to support label-free ultra-high throughput ambient MS and advance the state of the art for drug discovery.

Novel Aspect

Novel acoustic hyphenated sampling approach with mass spectrometry for ultra-high throughput label free screening and other applications in drug discovery.

References

1. Sinclair I. et al., J. Lab. Autom. 21, 19-26 (2016).
2. Van Berkel G., Kertesz V., Rapid Commun. Mass Spectrom. 29, 1749-1756 (2015).

MP-83 / EXPERIMENTAL MEASUREMENTS OF EFFECTIVE ELECTRIC FIELD RADII OF RECTILINEAR ION TRAPS

Chuanfan Ding - Fuxing Xu

Fudan University, Chemistry Department, Shanghai

Keywords: Linear ion trap, Electric field radius, Hyperbolic electrode, non-hyperbolic electrodes.

Introduction

The electric field radius is a fundamental character of an ion trap, it measures the basic properties of ion trap. For a quadrupole ion trap with idea hyperbolic electrode, its electric field radius can be simply obtained by its geometrical structure, but it is very difficult to calculate the electric field radius if it was built with non-hyperbolic electrodes. In this work, the electric field radii of different size rectilinear ion traps were measured by experiments.

Methods

The rectilinear ion traps with 5.0 mm x 5.0 mm, 5.50 mm x 5.00 mm, 6.00 mm x 5.00 mm, 6.50 mm x 5.00 mm, and 7.00 mm x 5.00 mm respectively were used for radius testing. All experiments were performed using an electrospray ionization(ESI) rectilinear ion trap mass spectrometry. The digital power supply was used for both ion trapping and dipolar excitation of mass selected ion during CID experiments. And the electric field radius of each ion trap was deduced by the results.

Results

For a quadrupole linear ion trap (LIT)with hyperbolic electrodes, the ion motion inside ion trap can be described by Mathieu equation.

$$q = \frac{4eV}{mr^2\Omega^2} \quad (1)$$

And we can get following (2) from (1):

$$\frac{m}{e} = \frac{V}{qr^2\pi^2} T^2 \quad (2)$$

Whether an ion of m/e can be trapped depends on the radial size r , or it is so called the electric field radius, the RF angular frequency $\Omega=2\pi f$, the voltage (V), and q . The property of an ion trap, such as the ion trapping capability, ion mass range and so on are mainly dependent on the ion trap parameters, And particularly, the electric field radius. For a LIT with non-hyperbolic electrode such as rectilinear ion trap, its r cannot be directly deduced from equation (1). But it can be measured by CID experiment when a ion of m/e was performed at the optimal RF period of T according to (2) when V is fixed.

For the rectilinear ion traps with 5.0 mm x 5.0 mm, 5.50 mm x 5.00 mm, 6.00 mm x 5.00 mm, 6.50 mm x 5.00 mm, and 7.00 mm x 5.00 mm, their measured electric field radii were measured, and they are 4.89 mm, 5.15 mm, 5.30 mm, 5.79 mm and 6.10 mm respectively.

Conclusions

The electric field radii of rectilinear ion trap can be measured by experiments. For proving the reliability of this method, a hyperbolic electrode ion trap with $r_0=4.00$ mm was built and its measured electrode field radius was 4.07 mm.

Novel Aspect

The electric field radii of different size rectilinear ion traps were measured by experiments.

MP-84 / UNIVERSAL ELECTRODE FOR QUADRUPOLE ELECTRODE SYSTEMS

Mikhail Dubkov - Mikhail Burobin - Vladimir Ivanov - Alexandr Malutin - Victor Gurov

Ryazan State Radio Engineering University, Department of General and Experimental Physics, Ryazan

Keywords: quadrupole mass filter, universal hyperbolic electrode

Introduction

It is offered to create quadrupole mass-analyzers as short identical units-monoblocks with the relevant electrodes concatenated mechanically and electrically to form the unified system of electrodes [1, 2]. Each unit-monoblock consists of all-purpose hyperbolic electrodes which are concatenated with the help of ceramic isolators.

Methods

To create an all-purpose hyperbolic electrode the electrolytic molding technology [3] which allows making thin-walled copper electrodes (1–2 mm) with high accuracy was applied.

Results

The design principle described is worked in practice. The all-purpose hyperbolic electrode is made of copper by the electrolytic molding method, its length is 75 mm, its wall thickness is 1–1,5 mm. The electrode is clothed in coating of 5–10 micron thickness by means of electrical deposit. The peripheral zone houses manufacturing holes which allow fastening electrodes together with a screw through ceramic isolators. The webbing along the electrodes' ends in addition to its intended purpose is for the fastening of separate electrodes' units together during the assembling of extended systems as well as for the setting of

the whole electrodes' system in vacuum chamber. The quadrupole mass filter was created with the using of the all-purpose hyperbolic electrodes.

Conclusions

The all-purpose hyperbolic electrodes allow creating quadrupole electrodes' system as separate blocks which can be connected with each other to get transit-time mass-analyzers of practically any length. The longitudinal distortion of electrical field doesn't occur as the accuracy depends on the performance quality of a separate monoblock.

Novel Aspect

The application of the all-purpose hyperbolic electrode created by the electrolytic molding method with one precision form usage for quadrupole mass-analyzers creation.

References

1. Gurov V. S., Dubkov M. V., Burobin M. A., Kharlanovl. A., Vestnikof RSREU, 52, 50-154 (2015).
2. Gurov V. S., Dubkov M. V., Burobin M. A., Ivanov V. V., Radiotekhnika, 8, 81-85 (2016).
3. Sheretov E. P., Gurov V. S., Dubkov M. V. e. a. Patent RF, no. 2091902, 1996.

MP-85 / DEGRADATION STUDIES OF THE ANTIBIOTIC SULFAMETHOXAZOLE UTILIZING HPLC/MS

Kevin Eckey - Uwe Karst

University of Münster, Institute of Inorganic and Analytical Chemistry, Münster

Keywords: Sulfamethoxazole, HPLC, ICP-MS, Hyphenation, Environmental

Introduction

The removal of micropollutants from surface waters is a topic of growing awareness. Not all xenobiotics are sufficiently removed by wastewater treatment plants, so additional disinfection steps such as ozonation or direct photolysis utilizing UV light are investigated. Many of these procedures are already well elucidated regarding degradation efficiency or energy costs but only few studies deal with the identification of formed products and their toxicity.

Methods

Sulfamethoxazole (SMX) solutions adjusted to pH 3 and pH 8 respectively were degraded using direct photolysis (220-500 nm). The resulting transformation product mixture was characterized using HPLC/ESI-HRMS and fragmentation experiments. Sulfur-containing compounds were quantified using HPLC/ICP-TQ-MS.

Results

In the present study, the degradation of the frequently used antibiotic sulfamethoxazole (SMX) was investigated using photolysis as degradation process focusing on the identification and characterization of formed transformation products by means of high performance liquid chromatography coupled to electrospray-mass spectrometry (HPLC/ESI-MS). The SMX solutions were pH-adjusted using a phosphate buffer to study the pH dependency of the product generation. In a further step, formed transformation products containing sulfur were quantified using HPLC/ICP-TQ-MS. From the quantitative data, rate constants of the photolytic degradation could be determined.

Conclusions

More than 30 degradation products were found, with some being already described in literature. It was shown that the pH value has a great impact on the SMX degradation behavior and on the selection of

formed transformation products.

Novel Aspect

An HPLC/ICP-TQ-MS hyphenation was developed for the quantification of sulfur-containing organic compounds and the method was optimized to an LOD of about 100 ng/mL.

MP-86 / THE IMPORTANCE OF ENHANCED MOLECULAR IONS IN GC-MS

Alexander B. Fialkov⁽¹⁾ - ***Uri Keshet***⁽¹⁾ - ***Tal Alon***⁽¹⁾ - ***Aviv Amirav***^(1,2)

School of Chemistry, Tel Aviv University, Tel Aviv (1) - Aviv Analytical Ltd, 3 Haarad Street, Tel Aviv 69710, Israel (2)

Keywords: Benefits of enhance molecular ions

Introduction

A major limitation of GC-MS is the weakness or absence of molecular ions in its standard EI mass spectra. Enhanced molecular ions by the use of EI of cold molecules (named "Cold EI") has many benefits.

Methods

We used GC-MS with Cold EI based on the combination of an Agilent 5977 MSD with the Aviv Analytical supersonic molecular beam (SMB) interface and fly-through ion source with electron ionization of internally cold molecules in the SMB (hence the name Cold EI) and the following are its benefits:

Preliminary data

Cold EI provides both enhanced molecular and fragment ions thus it is the only "soft" ionization method that is compatible with library identification.

Unexpectedly, enhanced molecular ion actually improves NIST library identification probabilities.

Molecular ions strongly decline with mass in standard EI (20% for added carbon in hydrocarbons) while in Cold EI it is about size independent.

Cold EI eliminates vacuum background and reduces column bleed thus make it easier to visually identify the molecular ions than standard EI and/or ESI-LC-MS.

Enhanced molecular ions and high mass fragments provide extended structural and isomer information.

Molecular ions enables the use of TAMI software with unit resolution quadrupole MS for the conversion of isotopomeric abundances into elemental formula.

Matrix interference is reduced with mass thus it is reduced on molecular ions no less than with MS-MS on fragment ions.

The molecular ions in large compounds provide superior sensitivity in SIM or RSIM.

Novel Aspect

How enhanced molecular ions in GC-MS with Cold EI improve GC-MS based sample identification and several other GC-MS performance aspects

MP-87 / GC-MS WITH COLD EI - EXTENDING THE RANGE OF COMPOUNDS AND APPLICATIONS AMENABLE FOR GC-MS ANALYSIS

Aviv Amirav - ***Alexander B. Fialkov*** - ***Tal Alon*** - ***Uri Keshet***

Tel Aviv University, School of Chemistry, Tel Aviv

Keywords: GC-MS with Cold EI

Introduction

The major shortcoming and Achilles Heel of GC-MS is its inability to analyze relatively non-volatile and thermally labile compounds. This problem is further exacerbated via the often weakness or absence of molecular ions in standard EI-MS.

GC-MS with Cold EI is based on GC and MS interface with a supersonic molecular beam (SMB) and electron ionization of the SMB vibrationally cold sample compounds in a fly-through ion source.

Methods

We found that GC elution temperatures can be significantly lowered upon the reduction of column length and increase of carrier gas flow rate. The increase of column flow rate also lowers the elution temperatures from the injector liner into the column while ion source degradation is inherently eliminated in the fly-through ion source. Thus, GC-MS with Cold EI provides significant extension of the range of compounds and applications amenable for analysis.

Preliminary data

The transition from simple to incompatible for GC-MS analysis. As the sample compound size and polarity is increased its analysis by standard EI gradually becomes more difficult in all the performance aspects: TIC signal is reduced, molecular ion abundance is reduced, baseline noise is increased, peak tailing is increased, ion source degradation is increased, the response becomes non-linear and analysis time is increased until the analysis becomes non-practical. In contrast, with Cold EI these performance aspects decline are much weaker.

Large and labile compounds analysis. The analysis of all organic explosives including TATP, PETN, Tetryl and HMX will be demonstrated plus thermally labile carbamate pesticides, underivatized steroids and large triglycerides in oil and blood plus large synthetic organic compounds.

Isomer distribution analysis for fuels and oils characterization. The availability of molecular ions to all hydrocarbons uniquely enables isomer distribution analysis for improved fuels and oils characterization.

Novel Aspect

How GC-MS with Cold EI extend the range of compounds amenable for GC-MS analysis with new demonstrated applications

MP-88 / STUDY OF ENDOCRINE DISRUPTING COMPOUND RELEASE FROM MEDICAL DEVICES IN PLASMA SAMPLES OF PREMATURE NEWBORNS THROUGH AN ON-LINE SPE UHPLC-MS/MS METHOD

Fabio Gosetti⁽¹⁾ - *Bianca Bolfi*⁽¹⁾ - *Marcello Manfredi*⁽¹⁾ - *Elisa Robotti*⁽¹⁾ - *Federica Ferrero*⁽²⁾ - *Gianni Bona*⁽³⁾ - *Emilio Marengo*⁽⁴⁾

University of Piemonte Orientale, Department of Science and Technological Innovation, Alessandria (1) - Ospedale Maggiore della Carità, Terapia Intensiva Neonatale, Alessandria (2) - University of Piemonte Orientale, Department of Health Sciences, Novara (3) - Alessandria, Department of Science and Technological Innovation, Alessandria (4)

Keywords: endocrine disrupting compound release; migration; medical devices; mass spectrometry; perfluorochemicals

Introduction

Endocrine disrupting compounds (EDCs) are exogenous substances that mimic/antagonize the endogenous hormone effect. Everyday people are exposed to many products containing EDCs such as plastic containers, cladding films in aluminum food cans, medical devices (MDs), and personal care products. The

underdeveloped reproductive apparatus and immune system of newborns can be exposed to unusually high concentrations of EDCs difficult to metabolize or excrete [1].

Methods

The study deals with the development of an online UHPLC-MS/MS method for the determination of 25 EDCs belonging to different chemical classes (perfluorochemicals PFCs, bisphenols, alkylphenols, phthalates, and parabens). To evaluate the possible EDC migration from different MDs of diverse materials used for blood sampling or for parenteral nutrition therapies, the devices were put in contact with physiological and parenteral solutions at prefixed time.

Results

The first step of the study was to identify a possible release of EDCs from the MDs and successively the EDC accumulation in plasma of the hospitalized newborns during the therapies to which infants have been subjected.

The results demonstrated the released of some EDCs (such as parabens at concentration of few pg mL^{-1} , diethylphthalate and dibutylphthalate at ng mL^{-1} level), when using physiological solution as the leaching agent. In addition, the presence of a PFOS structural isomer was found at 25.0 ng mL^{-1} : this compound was released in part from the infusion tubes and mostly from the venous catheters.

The releases from the infusion tubes of the same materials (PE) from different manufacturers are very comparable, whereas the catheter of silicone releases greater amount of EDCs than that of polyurethane.

As concerns plasma samples, butylparaben has greater concentrations in healthy newborns than in premature ones, while PFOA and PFHpA show greater concentrations in newborns hospitalized in neonatal intensive care unit.

Conclusions

The developed and validated method is attractive for routine analysis as it is characterized by high sensitivity, selectivity, robustness and reduced runtime analysis.

The release from MDs of different plastic materials were investigated and diverse amount of EDCs ranging from a few pg mL^{-1} to about twenty ug mL^{-1} were found. EDCs were present also in the plasma of newborns and their monitoring have to be taken into consideration to avoid pathologies.

Novel Aspect

A new method has been developed and validated to quantify the migration of 25 EDCs belonging to different chemical classes (including PFCs) from MDs in plasma.

References

1. Schug T. T., Janesick A., Blumberg B., Heindel J. J., Journal of Steroid Biochemistry & Molecular Biology, 127, 204 (2011)

MP-89 / TRIODE-TYPE MCP-BASED ION DETECTOR ENABLING MINIATURE MASS SPECTROMETERS TO OPERATE AT LOW VACUUM UP TO 1PA

Masahiro Hayashi - Yuya Washiyama - Kentaro Ishii - Joji Sakakiyama

Hamamatsu Photonics K. K., Electron Tube Division, Hamamatsu

Keywords: microchannel plates, triode structure, low vacuum operation, miniature mass spectrometer

Introduction

Recent interest in portable mass spectrometry applications has spurred a great deal of research into the development of miniature mass spectrometers (MS) such as desktop type explosives trace detectors

(ETDs). In the near future, there is a possibility that ETDs used at checkpoint screening will change from the current ones based on ion mobility spectrometers to ones based on next-generation MS.

In order to retain the major advantage of portability, miniature MS systems typically use small vacuum pumps with relatively low pumping speeds, meaning the miniature MS has to be operated at relatively low vacuum levels ($>1e-2$ Pa). Also in terms of sensitivity and vacuum condition, almost all miniature MS systems under development have adopted the ion trap analyzer. In these system configurations, ion detectors are required to have both low vacuum operation ($>1e-2$ Pa) and high gain ($>1e+6$). Currently, channel electron multipliers (CEM) are mainly used in miniature MS, but it is difficult for CEM to satisfy the required specifications.

Conventional electron multiplier based detectors including CEM and microchannel plates (MCP) cannot be adapted to operate at low vacuum ($>1e-2$ Pa) because of ion feedback (IFB), which causes discharges and decreases in S/N. However, the reduced dimensions of a typical MCP detector, which can be less than the mean free path in conventional electron multiplier based detectors, make it suitable for operation at low vacuum. The major source of IFB, which is a peculiar phenomenon of secondary electron amplification, was found to be the creation of ions from residual gas molecules.

We concluded that controlling the residual gas ions is more important than suppressing the ion generation. In order to realize this idea, we introduced a triode structure into the conventional MCP detector to enable it to operate at low vacuum condition.

Methods

The mechanism of IFB in an MCP detector was investigated, and it was found that 1) IFB is caused by EI-generated residual gas ions, which are generated between the MCP-out and an anode by the reaction of residual gas and electrons from the MCP, 2) after traveling in reverse through the MCP channel, the residual gas ion generates secondary electrons by impinging on the MCP channel surface, and 3) IFB occurs repeatedly like an avalanche phenomenon and it becomes a huge noise signal. We concluded that controlling residual gas ions is more important than suppressing ion generation.

Almost all MCP-based ion detectors have a bi-planar structure which consists of a chevron stack MCPs and an anode. The potential is MCP-in $<$ MCP-out $<$ anode, so the residual gas ions(+) return to the MCP and cause an IFB. The triode structure has a mesh electrode placed between MCP-out and a dynode as the third electrode. The aim of adopting the triode structure is to prevent positive ions returning to the MCP. Therefore, the novel potential is MCP-in, dynode $<$ MCP-out $<$ mesh, so that the dynode can function as an ion capture electrode and the mesh as an electron capture electrode.

It is confirmed that the triode structure prevents ion feedbacks at low vacuum condition by evaluating various characteristics. Since the usual electron gun could not be operated at low vacuum condition, ultraviolet rays were used as input signals instead of electrons. The vacuum pressure of the chamber which contained the triode MCP detector was controlled by opening and closing the manual gate valve and monitored with a vacuum gauge. It was difficult to keep the vacuum pressure constant in the high pressure region, and it slowly deteriorated during the measurement.

Results

The prototype triode MCP detector was examined using standard MCP detector characterization equipment in high vacuum levels ($1e-4$ Pa or less), which is the recommended environment for the standard MCP detector. The gain characteristic of the triode MCP detector was $5e+6$.

The triode detector operated in two potential modes: a) bi-planar mode [MCP-in $<$ MCP-out $<$ mesh = dynode (connected)], b) triode mode [MCP-in = dynode $<$ MCP-out $<$ mesh]. When the gain of MCPs was lower ($1e+3$ or less), IFB did not occur in both modes. When the gain was higher ($1e+5$ or more), IFB occurred at $1e-3$ Pa in bi-planar mode without input signals, but did not occur in triode mode up to 1 Pa. It was found that the MCP detector was capable of high pressure operation by adapting the triode structure with a novel potential mode.

Since the charge of the input signal could not be measured by using ultraviolet rays, the peak channel from PHD characteristics analysis using a multichannel analyzer system, which measured the amount of charge from the MCP output signal, indicates the gain at each degree of vacuum level. Despite more than $1e-2$ Pa,

the dark noise which derived from the IFB could not be confirmed in the PHD characteristics. When the vacuum pressure was 1 Pa, the calculated gain of $5e+6$ was maintained.

In order to confirm the stability of the detector operation, the outputs of the detector were monitored for about 90 minutes while input signal was turned on and off at 15-minute intervals. The output of the triode detector was stable from the long-term examination results irrespective of the degree of vacuum.

Conclusions

The MCP detector for low vacuum operation was developed by combining the triode structure with a novel potential mode. A compact detector for miniature MS could be realized by simplifying the structure in terms of supply voltage conditions and achieve the target gain of $1e+6$ at 0.1 Pa. We will confirm the MCP detector operation in a lower vacuum region by improving the current evaluation system, and optimize the MCP detector for a commercialized MS.

Novel Aspect

The MCP detector for low vacuum operation was developed by combining the triode structure with a novel potential mode. A compact detector for miniature MS could be achieved with the target gain of $1e+6$ at 1 Pa.

MP-90 / THE DYNAMIC BEHAVIOR OF NANO LIQUID BRIDGE DURING AMBIENT SAMPLING IONIZATION IN SCANNING PROBE ELECTROSPRAY IONIZATION MASS SPECTROMETRY

Bui Kamihoriuchi⁽¹⁾ - ***Yoichi Otsuka***⁽¹⁾ - ***Aya Takeuchi***⁽¹⁾ - ***Futoshi Iwata***⁽²⁾ - ***Takuya Matsumoto***⁽¹⁾

Osaka Univ., Department of chemistry, Toyonaka (1) - Shizuoka Univ., Department of engineering, Hamamatsu (2)

Keywords: Tapping-mode scanning probe electrospray ionization (t-SPESI), Nano liquid bridge, Mass spectrometry imaging, Atmospheric pressure ionization, Ambient sampling ionization (ASI)

Introduction

Tapping-mode scanning probe electrospray ionization (t-SPESI) have been used for mass spectrometry imaging. In t-SPESI, the formation of a liquid bridge and the generation of an electrospray by a vibrating capillary probe realizes the ambient sampling ionization (ASI) [1-4]. However, the dynamic behavior of the charged liquid bridge at the probe apex is unclear. In the present study, we studied a probe motion during ASI with the optical detection method.

Methods

We developed the measurement system of the probe displacement. The optical lever method and the sample stage controller was incorporated in the experimental setup of t-SPESI. Mass spectrometry was conducted at the same time to capture the time changes of ion intensity. The commercial available silica emitter and the Rhodamine B thin film on glass was used as a probe and the specimen. The physical properties of the probe were simulated with finite element method.

Results

The spring constant of the probe was estimated to be about 37 N/m (18 mm long probe) and 5.0 N/m (35 mm long probe). We found that the probe with lower spring constant is suitable for the measurement of the probe displacement during ASI with high sensitivity. The probe apex was drawn toward the sample just before the contact, indicates the timing of liquid bridge formation. During the retraction of probe from the sample, the probe was released from a liquid bridge. Considering the retraction distance and the simulated spring constant, the force applied to the apex of the probe was about 25 μ N just before the breakage of liquid bridge. At that time, the ion intensity became stronger and it gently became weak within 1 to 2

second. From these result, it was found that the sampled molecules were ionized over time just after the breakage of liquid bridge. In addition, the result of experiment with a nanopipette (a nanoscale diameter opening glass pipette) will also be presented.

Conclusions

New optical lever method was developed for SPESI. The formation and the breakage of the liquid bridge was captured by the displacement of probe during ASI. The sampled molecules are ionized over several seconds after the breakage of liquid bridge. The pulling force by the probe to break the liquid bridge was estimated to be larger than the expected surface tension of liquid bridge.

Novel Aspect

This report is the first result to capture the dynamic probe motion during the sampling and ionization process using liquid bridge and electrospray in SPESI.

References

1. Y. Otsuka et. al., Rapid Commun. Mass Spectrom. 26, 2725 (2012).
2. Y. Otsuka et. al., Analyst 139, 2336 (2014).
3. Y. Otsuka et. al., J. Mass Spectrom. 50, 1157 (2015).
4. T. Kohigashi et. al., Mass Spectrometry (Tokyo) 5, S0054 (2016).

MP-91 / INTEGRATION OF A DROPLET MICROFLUIDIC CHIP WITH MASS SPECTROMETRY – A STEP TOWARDS HIGH THROUGHPUT CELLULAR ANALYSIS

Emily Kempa⁽¹⁾ - ***Xin Li***⁽²⁾ - ***Clive A. Smith***⁽²⁾ - ***Perdita Barran***⁽¹⁾

University of Manchester, Manchester Institute of Biotechnology, Manchester (1) - Sphere Fluidics Limited, Chemistry, Cambridge (2)

Keywords: Droplets, Microfluidics, lab-on-a-chip, Cells

Introduction

A major bottleneck to directed evolution workflows is the throughput required to screen a large number of genomic variants for the most efficient biocatalysts [1]. Current screening methods often rely upon well-plate fluorescent readouts which become problematic when the product of interest is not inherently fluorescent; hence, 'label-free' detection is highly desirable.

Methods

Mass Spectrometry coupled with a droplet microfluidic chip system offers a 'label-free' screening platform in which each genomic variant of the combinatorial library is encapsulated in a ~500 picolitre droplet suspended within an immiscible fluoruous oil phase. MS analysis of the cell filled droplets allows for the identification and capture of the highest product producing genomic variants.

Results

Droplet generation occurs at a rate of 5 Hz within a custom designed microfluidic chip followed by immediate flow of the droplets into a mass spectrometry source for ionization and then subsequent analysis. Loss of droplet order and droplet merging is avoided through 'plugs' of the oil phase and hence, each droplet arises in the TIC as a distinct peak at the same rate as generation. This allows for a mass spectrum for each peak to be obtained such that information regarding products or metabolome fingerprint in 1 discreet droplet can be extracted. This has been performed for a range of standards (including droplet mixtures [2]) and for bacterial cells cultured within droplets as proof of concept. Other on-chip processes such as droplet splitting have also been successfully trialed independently, with the aim to incorporate droplet capture of those identified by the MS as 'hits'.

Conclusions

We demonstrate direct coupling between chip and two different Ion-Mobility enabled Mass Spectrometers to analyse discreet droplets without the prior removal of the contents (both protein and small molecule standards and cells) from the oil phase. This has been achieved with a throughput of ~5 Hz, which already surpasses the rate of many current sample introduction systems.

Novel Aspect

Novel aspects include ultra-high throughput capabilities (~5 Hz) and miniaturization of sample. Introduction platforms for life science applications.

References

1. Aharoni A., Griffiths A. D., Tawfik D. S., *Curr. Opin. Chem. Biol.*, 9, 210-216 (2005).
2. Smith C. A., Li X., Mize T. H., Sharpe T. D., Graziani E. I., Abell C., Huck W. T. S, *Anal. Chem.*, 85, 3812-3816 (2013).

MP-92 / PORTABLE GAS CHROMATOGRAPHY MASS SPECTROMETER WITH A RAPID HEATING CARBON NANOTUBE PASTE

Bongyoon Yi⁽¹⁾ - *Seung Yong Kim*⁽²⁾ - *Wanseop Jeong*⁽¹⁾ - *Jae Yeong Eo*⁽²⁾ - *Byoung Seob Lee*⁽²⁾ - *Byeongwon Kang*⁽³⁾ - *Hyun Sik Kim*⁽²⁾

Korea Basic Science Institute/Mass Spectrometry & Advanced Instrumentation Group, Chungbuk National University/Department of Physics, Cheongju (1) - Korea Basic Science Institute, Mass Spectrometry & Advanced Instrumentation Group, Cheongju (2) - Chungbuk National University, Department of Physics, Cheongju (3)

Keywords: portable GCMS, cold electron, nanovalve, carbon nanotube, VOC

Introduction

Recently, a miniaturized ion trap mass spectrometer with cold electron ionization was developed by our group as an ultra-portable mass spectrometer [1]. To extend the applicable field of the portable MS to a complex gas mixture analysis, we integrated a gas chromatography (GC) system to the miniaturized ion trap mass spectrometer using a nano-liter valve interface and evaluated the analytical performances of the system using typical volatile organic compounds.

Methods

A new portable GC was developed using a new carbon nanotube (CNT) heating paste. The GC column was heated by a bobbin type heater made of anodized aluminum and CNT heating paste applied on the inner surface of the heater. The miniature cylindrical ion trap MS was developed using a cold electron source to ionize analyte gas. An accurate nano valve was developed and used as an interface between GC and the MS to avoid the failure of ion pump by abrupt increase of chamber pressure.

Results

The analytical characteristics of the portable GCMS were evaluated as individual components and as complete GCMS system. The performances of gas chromatography were evaluated with an adsorption concentrator using standard VOC samples while changing the column temperature gradient up to 100 degree/sec for rapid analysis. The performances of cylindrical ion trap were tested with nanoliter injection valve using the standard gas. With built in 2 L/sec ion pump and nanoliter injection valve, maximum 10 times/min sample injection was the limiting factor of analysis frequency while maintaining the operable pressure of the ion trap. The electron ionization was using electrons from electron multiplier induced by UV photons from a diode. The optimal transmission condition of the electrons to the ion generation region was

derived from SIMION simulation and used for the analytical performance test of ion trap. The developed GC/MS system was less than 10 kg and showed better than unit MS resolution and sub-ppm detection limits for typical VOCs.

Conclusions

Using a portable gas chromatography and a miniaturized cylindrical ion trap mass spectrometer, a new portable GCMS was introduced with CNT heating paste showed a high temperature-gradient for rapid analysis, and with a cold electron ionization source optimized for higher electron transmission to the ion generation region showed improved detection efficiencies. The integrated system showed detection limit of down to sub-ppm levels for the standard VOCs.

Novel Aspect

A novel hand-carriable GCMS was developed using a carbon nanotube heating paste for GC and cold electron ionization for portable ion trap MS to detect VOCs realtime in the field.

References

1. H. S. Kim, S. Y. Kim, and M. Yang, *Anal. Chem.* 2012, 84, 3635-3639

MP-93 / ANALYTICAL PROCEDURE FOR DETERMINATION OF PERFLUOROALKYL ACIDS, PARABENS AND COTININE FROM A LOW VOLUME HUMAN SERUM SAMPLE BY LC-MS/MS

Jani Koponen - Hannu Kiviranta

National Institute for Health and Welfare (THL), Health Security, Kuopio

Keywords: biomonitoring, cotinine, epidemiological study, paraben, perfluorinated substances

Introduction

In studies on the health effects of environmental chemicals it is desirable that multiple compound groups can be analyzed from the same low volume serum aliquot. Thus, we have developed and validated a single analytical procedure for determination of three different compound groups, perfluorinated substances (PFAS), parabens and cotinine (as a metabolite of nicotine), which are associated with a variety of adverse health outcomes [1-7].

Methods

Serum sample of 0.1 mL was pre-treated with two step extractions, and all the analytes were separated and detected with LC-MS/MS. The limit of quantitation (LOQs) for PFAS, parabens and cotinine were 0.1 - 0.3 ng/mL, 0.20 - 1.0 ng/mL and 0.10 ng/mL, respectively. These limits are highly sufficient for studies where the exposure effects of PFASs, parabens and tobacco smoke is evaluated.

Results

Validation data shows that the method has excellent accuracies and repeatabilities. For PFAS the accuracies ranged from 94 to 106% and the repeatabilities (as RSD%) from 2.2 to 5.6 %. The accuracies and RSDs for parabens were 73 - 120 % and 2.2 - 9.7 %, respectively, and for cotinine 102 - 106 % and 2.9 - 3.5 %. The main PFAS found in population based serum were perfluorooctanesulfonate (PFOS, range 1.3 - 23 ng/mL), perfluorooctanoate (PFOA, range 0.54 - 6.6 ng/mL), perfluorononanoate (PFNA, range 0.19 - 4.1 ng/mL) and perfluorohexanesulfonate (PFHxS < 0.15 - 5.5 ng/mL). For parabens the methyl (range < 1.0 - 49 ng/mL) and propyl (range < 0.2 - 1100 ng/mL) conjugates were the abundant compounds found in this study. Cotinine level in serum of tobacco smokers was over 50 ng/mL, whereas in serum of non-smokers the level was usually below 0.1 ng/mL.

Conclusions

A simple, sensitive and high throughput method for analysis of 21 PFAS, 8 parabens and cotinine from 0.1 mL of serum was developed and validated. Analysis throughput time for 150 serum samples is one week. Besides the high throughput performance, the method is suitable for studies where the sample volume is limited since only a low sample volume is needed.

Novel Aspect

This is the first method to gather these 3 classes of environmental pollutants (perfluorinated substances, parabens, and cotinine) in a single procedure.

References

1. Bonfeld-Jorgensen EC., Long M., Bossi R., Ayotte P., et al., *Environ Health*, 10, 88(2011)
2. Nelson JW., Hatch EE., Webster TF., *Environ Health Perspect*, 118, 197(2009)
3. Steenland K., Tinker S., Frisbee S., Ducatman A., Vaccarino V., *Am J Epidemiol*, 170, 1268 (2009)
4. Frisbee SJ., Shankar A., Knox SS., Steenland K., et al., *ArchPediatrAdolescent Med*, 164, 860 (2010)
5. Knox SS., Jackson T., Frisbee SJ., Javins B., Ducatman AM., *J ToxicolSci*, 36, 403 (2011)
6. Giulivo M., Lopez de Alda M., Capri E., Barceló D., *Environ Res*, 151, 251(2016)
7. Johnson GK., Squier CA., *Pediatr Dent*, 15, 169 (1993)

MP-94 / COLD VAPORIZATION OF TISSUES BY PICOSECOND INFRARED LASER ABLATION (PIRL) - UNIQUE ACCESS TO THE ORIGINAL PROTEOFORM COMPOSITION.

Marcel Kwiatkowski⁽¹⁾ - ***Hartmut Schlüter***⁽²⁾

University of Groningen, Department of Pharmacokinetics, Toxicology and Targeting, Groningen (1) - University Medical Center Hamburg-Eppendorf, Institute of Clinical Chemistry, Mass Spectrometric Proteomics, Hamburg (2)

Keywords: Picosecond infrared laser ablation (PIRL); Protein extraction; Tissue proteomics; Proteoform; Clinical proteomics

Introduction

The proteoform composition of tissues has a high diagnostic value. Proteoforms regulate biological processes and an aberrant composition may cause diseases. Intact proteoforms are exposed to enzymatic degradation reactions during sample preparation and methods are needed that provide access to the original proteoform composition of tissues. PIRL ablation minimizes degradation reactions and provides access to the original proteoform composition of tissues.

Methods

PIRL-ablated proteins were analyzed by LC-ESI-MS/M and MALDI-MS to investigate if their chemical composition and enzymatic activity remained unaltered. Tissues were homogenized either by conventional mechanical methods or by PIRL ablation. The homogenates were separated by 2-DE and SDS-PAGE, followed by quantitative LC-MS/MS analysis. Enzymatic degradation reactions were analyzed by proteoform SDS-PAGE migration profiles.

Results

Comprehensive MS analysis showed that the exact chemical composition of proteins was not altered during the PIRL ablation process and even posttranslational modifications and enzymatic activities were reserved (1). 2-DE analysis indicated that PIRL yielded a higher number of intact proteoforms compared to mechanical homogenization. These results were confirmed by SDS-PAGE migration profiles. Relative

degrees of proteolysis were significantly lower in case of PIRL (2%) compared to mechanical homogenization (23%) (2).

These results were further confirmed by spiking experiments. Total yield of intact alpha-casein proteoforms and number of identified phosphopeptides were significantly higher using PIRL, whereas mechanical homogenization showed significantly higher amounts of enzymatically degraded proteoforms.

PIRL-derived extracts are almost free of particles minimizing sample, thus providing a higher reproducibility and significantly higher protein identification rates compared to conventional tissue homogenization.

Conclusions

PIRL releases intact proteoforms from tissues without changing their exact chemical composition including posttranslational modifications. The ultrafast release and transfer of intact proteoforms into frozen condensates minimizes exposure to enzymatic degradation reactions and provides an unique access to the original proteoform composition of tissues.

Novel Aspect

PIRL provides an unique access to the original proteoform composition of tissues and facilitates molecular diagnostics of tissues at the proteoform level.

References

1. Kwiatkowski et al. *Angew. Chem. Int. Ed. Engl* 54(1): 285-288 (2015).
2. Kwiatkowski et al. *J Proteomics* 134: 193-202 (2016).

MP-95 / FROM NICHE APPLICATION TO ESSENTIAL RESEARCH TOOL: WHY SFC-MS IS THE MUST HAVE TECHNIQUE

G. John Langley - Julie Herniman

University of Southampton, Chemistry, Southampton

Keywords: Supercritical fluid chromatography - mass spectrometry

Introduction

SFC now delivers the chromatographic promises with robust and reliable instrument. The solvation power, selectivity and peak capacity of modern Ultrahigh Performance Supercritical Fluid Chromatography (UHPSFC) coupled with mass spectrometry affords the ideal analytical platform to address classes of compounds that were previously a challenge for GC-MS and RP-UHPLC-MS approaches.

UHPSFC-MS sits between LC-MS and GC-MS affording extended capability and complementarity and delivers unique solutions across a broad range of application areas, *e.g.* synthetic organics, pharmaceuticals, petrochemistry, lipids, nucleotides and many other areas.

Methods

Ultra-high performance supercritical fluid chromatography (UHPSFC) methods have been developed using a simple column screening protocol. A Waters Acquity UPC2 is coupled to a TQD mass spectrometer and a range of cosolvents and additives are utilised to deliver generic and bespoke methods. Further tuning of the organic make-up solvent has been used to deliver specific assays.

Results

UHPSFC, in the form of Ultra Performance Convergence Chromatography (UPC2) has been incorporated into the open access (walk-up and use) chromatography mass spectrometry facility within Chemistry at the University of Southampton in 2013. This technology complements existing liquid chromatography and gas

chromatography options (HPLC-MS (ESI/APCI), GC-MS (EI/CI) and affords high throughput reaction monitoring for new chemistries and highly specific research solutions, *e.g.* lipids, isocyanate polymers, petrochemistry, molecular wires and supramolecular chemistry.

Novel Aspect

New chromatography solutions, extension of open access mass spectrometry.

MP-96 / CONVENTIONAL-FLOW LIQUID CHROMATOGRAPHY–MASS SPECTROMETRY FOR EXPLORATORY PROTEOMIC ANALYSES

Juraj Lenčo

Department of Analytical Chemistry, Faculty of Pharmacy in Hradec Kralove, Charles University, Hradec Králové

Keywords: Conventional-flow LC-MS, exploratory proteomics.

Introduction

Although conventional-flow LC-MS platforms have penetrated targeted proteomics successfully, their possibilities in discovery-oriented proteomics have not yet been thoroughly explored. Our objective was to determine what optimization and adjustments to a conventional-flow LC-MS system must be undertaken to identify a comparable number of proteins as can be identified on a nanoLC-MS system.

Methods

A trypsin digest of HeLa cells lysate was used to compare proteomic identifications between nanoLC and conventional-flow LC-MS setup. Following parameters of the conventional-flow LC-MS configuration were optimized in order to improve identifications: DDA settings, column diameter, extra-column peak dispersion, temperature, flow rate, column length, mobile phase acidifiers and additives.

Results

We demonstrate that the amount of a complex tryptic digest needed for comparable proteome coverage can be roughly 5-fold greater, providing the column dimensions are properly chosen, extra-column peak dispersion is minimized, column temperature and flow rate are set to levels appropriate for peptide separation, and the composition of mobile phases is fine-tuned. Indeed, we identified 2 835 proteins from 2 µg of HeLa cells tryptic digest separated during a 60 min gradient at 68 µL/min on a 1 mm ×250 mm column held at 55°C and using an aqua–acetonitrile mobile phases containing 0.1% formic acid, 0.4% acetic acid, and 3% dimethyl sulfoxide.

Conclusions

Although we have not exploited all options for improving LC–MS performance, we conclusively show that conventional-flow exploratory proteomic analyses are feasible. Providing the instrumentation and method are adjusted appropriately, the extra sample amount need not correspond to the theoretical 178-fold greater burden [1].

Novel Aspect

Our results document that nanospray-free, conventional-flow LC–MS is an attractive alternative for bottom-up exploratory proteomics.

References

1. LenčoJ. et al., Analytical Chemistry, 90, 5381(2018).

MP-97 / MODERN MOBILE MASS SPECTROMETRY FOR ON-SITE ANALYSIS OF VOCs***Andrea Marcillo***⁽¹⁾ - ***Anja Widdig***⁽²⁾ - ***Claudia Birkemeyer***⁽¹⁾

University of Leipzig, Institute of Analytical Chemistry, Research Group of Mass Spectrometry, Leipzig (1) -
 University of Leipzig, Institute of Biology, Research Group of Behavioral Ecology, Leipzig (2)

Keywords: Mobile, GCMS, performance, VOCs, profiling**Introduction**

The fast development and exhaustive search for effective VOCs analysis with mobile MS has become an important research field, due to the advantages of miniaturized systems for on-site applications over laboratory-scale instrumentation. Nevertheless, it is important to rely on valid analytical platforms capable of achieving high sensitivity, reducing process times, and, specifically for on-site analysis [1], enabling multi-component determinations.

Methods

An extensive evaluation of three commercially available mobile GCMS, the Hapsite ER (Inficon, Germany) and the E2M(Bruker, Germany) with quadrupole mass analyzers, and the Torion T-9 (PerkinElmer, Germany) with a toroidal ion trap mass analyzer, was performed using a complex mixture of 24 compounds introduced by different devices, and compared to a state-of-the-art bench-top GC-2010 Plus coupled to a TQ-8040 triple quadrupole mass spectrometer (Shimadzu, Japan).

Results

First, the performance was assessed in terms of sensitivity, approximately only 50% of the analytes were detected in the mobile mass spectrometers independent of the introduction devices (thermal desorption TD, direct probe and SPME) and compared to the TD stationary system. In contrast, after average normalization, the reproducibility of consecutive measurements (n=3) seemed to be strongly influenced by the pre-concentration system. While ~30% RSD was determined for ab/adsorption-desorption methods (SPME and TD), even higher variance and unstable results were obtained using the direct probe system; all of them were considerably higher than the 15% RSD found on laboratory scale. In addition, spectra reproducibility was probed for its impact on spectral similarities with common reference spectra libraries, but still some of the corresponding softwares showed several limitations for data transformation, data transfer, and even selection of quantitation and confirming ions for proper determinations and statistical analysis.

Conclusions

Mobile devices for on-site applications of unknown samples still need to overcome several analytical limitations during qualitative and quantitative analysis, specifically for volatile profiling with multi-component determinations.

Novel Aspect

Mobile GCMS was assessed for appropriate method development and validation for analysis of complex standard mixtures and compared to state-of-the-art bench-top GCMS.

References

1. Kotiaho, T., J. Mass Spectrom., 31, 1–15 (1996).

MP-98 / DETERMINATION OF MULTI-RESIDUE ANTIBIOTICS IN SURFACE WATER BY SPE AND UHPLC-MS/MS

Derek J. Mattern⁽¹⁾ - ***Yiwen Li***⁽²⁾ - ***Zhiwei Gan***⁽²⁾ - ***Ruoying Liao***⁽²⁾ - ***Xiang Li***⁽³⁾ - ***Xiangdong Zhou***⁽³⁾ - ***Chengyuan Cai***⁽³⁾ - ***Yongming Xie***⁽³⁾ - ***Feng Qin***⁽⁴⁾ - ***Jingcun Wu***⁽⁴⁾

PerkinElmer, Analytical Science, Rodgau (1) - Sichuan University, Department of Environmental Science and Engineering, Chengdu (2) - PerkinElmer, Chromatography, Shanghai (3) - PerkinElmer, Chromatography, Woodbridge (4)

Keywords: UHPLC-MS/MS, antibiotics, solid-phase extraction, surface water

Introduction

Antibiotics are used as human and veterinary medicines to treat diseases. Studies have found various antibiotics in surface water, which has raised concern on the human health effects. To investigate the occurrence and fate of these low level contaminants and measure them in a variety of water matrices; sensitive, selective and reliable methods are needed. In this study, a selective solid phase extraction (SPE) method was evaluated for sample clean-up and analyte concentration for analysis on an UHPLC-MS/MS system.

Methods

Water samples were subjected to solid-phase extractions. Analytes were determined by coupling the UHPLC system to a highly sensitive PerkinElmer QSight 220 triple quadrupole mass spectrometer with dual ESI/APCI sources. Under the optimized conditions, the method was validated in terms of its specificity, sensitivity, linearity, accuracy, and precision. Recoveries from sample matrix were evaluated by fortifications of analytes at 3 different concentrations.

Results and Conclusions

In this study, to effectively remove interfering components and enriched analyte of interest from complex environmental water sample matrices, several SPE extraction phases with different functionalities were evaluated, such as reversed phase, ion exchange and mixed mode cartridges. The extraction and elution conditions were also optimized for SPE procedures. Calibration curves were built with matrix-matched standards to compensate for any matrix effects. The regression coefficients (R²) for the calibration curves obtained for all antibiotics studied were all greater than 0.990, and the limits of quantification (LOQs) ranged from 0.1 to 5 µg/L for the targeted compounds. The mean recoveries from the samples spiked at three different concentrations, respectively, were between 70% and 120%. The intra-day and inter-day variations, expressed as RSD, were less than 5%, respectively. The results demonstrated that the accuracy and precision of the present method were acceptable for routine monitoring purposes. The analyte confirmation in sample was achieved by comparing their respective retention times and peak area ratios of two MRM transitions with those obtained from the corresponding reference standards. The method developed in this study provides an effective tool for determination of multi-antibiotics in surface water and a large number of actual water samples were analyzed using this method to evaluate the occurrence of the analytes in local water systems in Chengdu, China.

Novel Aspect

Highly selective SPE extraction method coupled to highly sensitive UHPLC-MS/MS for low level antibiotics analysis in surface water.

MP-99 / EVALUATION OF CAPILLARY ELECTROPHORESIS AND ION MOBILITY FOR CYCLIC LIPOPEPTIDES ANALYSIS

Andréa Mc Cann⁽¹⁾ - *Johann Far*⁽¹⁾ - *Cédric Delvaux*⁽¹⁾ - *Jaanus Liigand*⁽²⁾ - *Loïc Quinton*⁽¹⁾ - *Edwin De Pauw*⁽¹⁾

MSlab, Université de Liège, Liège (1) - Institute of Chemistry, University of Tartu, Tartu (2)

Keywords: Ion mobility - Capillary-electrophoresis – Lipopeptides – High resolution Mass Spectrometry – Reverse phase liquid chromatography

Introduction

Cyclic lipopeptides (CLPs) are produced by bacteria in a ribosome independent manner and consists in a hydrophilic peptide chain linked to a β -hydroxy fatty acid chain. CLPs are a very large and diverse group of molecules, that has recently gained lots of attention for their potential applications in agricultural or pharmaceutical industries [1, 2]. In that aim, the development of reliable analytical methods for new CLPs isolation and characterization is essential.

Methods

Capillary Electrophoresis-Mass Spectrometry was performed on a CESI 8000 instrument coupled to a tandem linear ion trap and high resolution mass spectrometry FT-ICR. A bare fused silica capillary of 90cm x 375 μ m (OD) x 50 μ m (ID) was used for the separation. In addition, CLPs were analyzed by ion mobility mass spectrometry on a timsTOF instrument.

Results

The obtained separation of up to four isobaric CLPs by capillary electrophoresis mass spectrometry will be presented. The fragmentation spectra enabled the identification of each compound detected. The mobility of the CLPs in solution will then be compared to the measured collision cross section in the gas phase by ion mobility mass spectrometry.

Conclusions

Capillary electrophoresis appears to be a reliable technique for cyclic lipopeptides separation and shows promising applications for new CLPs identification. Capillary electrophoresis and ion mobility mass spectrometry are complementary techniques for CLPs structural analysis.

Novel Aspect

To our knowledge, this is the first report on cyclic lipopeptides separation by capillary zone electrophoresis method.

References

1. Biniarz P., Łukaszewicz M., Janek T., Critical Reviews in Biotechnology, 37, 393-410 (2017)
2. Roongsawang N., Washio K., Morikawa M., International journal of molecular sciences, 12, 141- 172 (2011)

MP-100 / TANDEM IONISATION FOR GC(xGC)-MS: COMPLEMENTARY SOFT ELECTRON IONISATION FOR INCREASED CONFIDENCE IN COMPOUND IDENTIFICATION

Nick Bukowski⁽¹⁾ - ***Laura McGregor***⁽¹⁾ - ***Joe Blanch***⁽²⁾ - ***Pete Grosshans***⁽²⁾

SepSolve Analytical, Peterborough (1) - Markes International, Llantrisant (2)

Keywords: soft ionisation; TOF MS, tandem ionisation

Introduction

Electron ionisation (EI) at 70 eV is standard practice for GC-MS, due to the high degree of ionisation (and thus good sensitivity) and consistent mass spectra generated.

Nevertheless, the identification of individual compounds in complex samples remains challenging when multiple compounds in a chemical class have similar spectra, or weak molecular ions at 70 eV EI. This problem can be addressed by the use of soft ionisation to reduce the degree of ion fragmentation, but this

approach has been cumbersome to implement until now.

Methods

The mass spectrometer used in this study incorporates a unique ion source design enabling both hard and soft electron ionisation spectra to be acquired simultaneously – in a technique deemed Tandem Ionisation. Multiplexing between ionisation energies allows conventional 70 eV spectra to be obtained - for routine identification against commercial libraries - as well as simplified soft EI spectra for enhanced confidence in identification of compounds which may exhibit weak molecular ions and/or similar spectra at 70 eV.

Results

The soft EI spectra consist of only structurally-significant fragments for reduced demand on dynamic range and an increase in peak capacity. In addition, greater orthogonality between the mass spectra of isomeric compounds is evident, thus simplifying compound identification and reducing reliance on retention indices.

Conclusions

Fast multiplexing during acquisition greatly improves the performance of a TOFMS, enabling two complementary MS datasets to be acquired simultaneously.

Novel Aspect

The use of a novel, patented ion source design to provide both hard and complementary soft electron ionisation in a single run.

MP-101 / FLUOROUS MODIFICATION OF PEPTIDES FOR MASS SPECTROMETRIC APPLICATIONS

Adrienn Molnár⁽¹⁾ - ***Máté Berta***⁽²⁾ - ***Arnold Steckel***⁽³⁾ - ***Katalin Uray***⁽⁴⁾ - ***József Rábai***⁽²⁾ - ***Gitta Schlosser***⁽¹⁾

Department of Analytical Chemistry, ELTE Eötvös Loránd University, Budapest (1) - Department of Organic Chemistry, ELTE Eötvös Loránd University, Budapest (2) - Doctoral School of Chemistry, ELTE Eötvös Loránd University, Budapest (3) - MTA-ELTE Research Group of Peptide Chemistry, Hungarian Academy of Sciences, ELTE Eötvös Loránd University, Budapest (4)

Keywords: fluorous modification, chemical derivatization, LC-MS, polar peptides, proteomics

Introduction

Arginine- and lysine rich proteins are challenging targets in proteomics, due to the formation of short peptides during tryptic digestion. Small, polar peptides have very low retention in RP-LC-MS and are usually hidden under the standard conditions. Our aim was the fluorous derivatization of polar peptides to modify their chemical properties and to improve their chromatographic retention.

Methods

Mass spectrometric experiments were performed using a Bruker Esquire 3000+ ion trap mass spectrometer with Jasco PU-2085Plus HPLC. Peptides were reacted with perfluoroalkancarboxylic anhydrides $[F((CF_2)_nCO)_2O]$; where $n = 1 - 3$] in dry acetonitrile. The reaction was quenched with large excess of water. Products were analyzed without further purification.

Results

Model peptides containing various amino acids were synthesized by solid phase peptide synthesis: H-GZGG-OH and H-AZAA-OH [where Z = Gly/Ala, Arg, Cit, Ser]. Peptides were reacted with trifluoroacetic anhydride, pentafluoropropionic anhydride and heptafluorobutanoic anhydride, respectively. Propionic anhydride was used as reference. Under optimized conditions, conversion of the peptides to their *N*-acylated fluorous derivatives was complete in minutes. Effect of the various perfluoroalkyl groups was studied by reversed

phase LC-MS on a C18 column. MS/MS fragmentation patterns using collision-induced dissociation (CID) were also compared.

Conclusions

The terminal and side-chain amino group acylation by fluoros anhydrides is a very fast reaction under optimized conditions. Chromatographic properties of modified polar peptides were significantly improved. Beside the formation of sequential fragments, neutral loss of the fluoros acid was proved to be a characteristic fragmentation pathway.

Novel Aspect

Fluorous derivatization of polar peptides by anhydrides is a novel approach for improved detection of highly polar peptides by mass spectrometry.

Acknowledgement

Gitta Schlosser acknowledges the support of the MTA Premium Post-Doctorate Research Program of the Hungarian Academy of Sciences (HAS, MTA).

MP-102 / EVALUATING THE CAPABILITY OF DIFFERENT NON-TARGETED MS-BASED APPROACHES FOR THE COMPREHENSIVE CHARACTERIZATION OF OLIVE TREE-DERIVED MATRICES

Lucía Olmo-García⁽¹⁾ - *Nikolas Kessler*⁽²⁾ - *Aadil Bajoub*⁽³⁾ - *Carmen María Sánchez-Arévalo*⁽¹⁾ - *Arsenio Muñoz de la Peña*⁽⁴⁾ - *Heiko Neuweger*⁽²⁾ - *Carsten Baessmann*⁽²⁾ - *Alegría Carrasco-Pancorbo*⁽¹⁾ - *Alberto Fernández-Gutiérrez*⁽¹⁾

University of Granada, Department of Analytical Chemistry, Granada (1) - Bruker Daltonik GmbH, Bruker Daltonik GmbH, Bremen (2) - National School of Agriculture, Department of Basic Sciences, Meknès (3) - University of Extremadura, Department of Analytical Chemistry, Badajoz (4)

Keywords: non-targeted approaches; food metabolomics; LC-MS; GC-MS; complementarity platforms

Introduction

Virgin olive oil consumption has been associated to a lower incidence of some chronic diseases, largely due to its non-glyceridic minor compounds, which are also found in other vegetal tissues derived from the olive tree. Thus, the deep characterization of the main parts of the olive fruit and tree might be a key step when trying to find natural sources of bioactive compounds. To do it so, different MS-platforms were used (LC-ESI/APCI-MS and GC-APCI-MS).

Methods

Picudo cv. samples were prepared using an unselective liquid/solid-liquid extraction protocol. The extracts were analyzed by LC and GC coupled to a QTOF detector by means of ESI and APCI interfaces in the case of LC and an APCI source in GC. In LC, the analytes were eluted in a C18 column (2.1x 100 mm, 1.8 µm), with acidified water and acetonitrile at 40°C. The silylated extracts were also analyzed in GC (BR-5 column) with a T gradient from 150 to 320°C (ramped at 4°C/min).

Results

The extracts (from lyophilized olive pulp, olive seed, fruit skin, leaves and wood from the branches, as well as virgin olive oil, olive oil obtained from pitted and dehydrated fruits and oil obtained from the seed contained inside the pit) were prepared and analyzed. Data acquired were processed with MetaboScape® 3.0 (Bruker), which automatically extracts and combines isotopes, adducts and fragments belonging to the same compound into one feature. Information obtained from Metaboscape was studied to point out the adequacy of each methodology for the determination of each family of compounds, as well as to describe

the distribution of the identified plant metabolites within the selected olive oils and tissues (their presence, absence, relative area in each matrix and relative response in each platform were checked). The identification of more than 130 compounds from the extracts (phenolic compounds, triterpenic acids and dialcohols, tocopherols, sterols and free fatty acids) was achieved.

Conclusions

- The potential of different non-targeted MS-based approaches to cover the metabolome of 8 *Olea europaea* derived matrices was assessed.
- The coupling to a high resolution MS through different ionization sources (ESI and APCI in LC and APCI in GC) allowed the detection of compounds from different chemical classes.
- The use of a great number of standards and the annotation strategies (MetaboScape) lead to the identification of more than 130 compounds in the extracts.

Novel Aspect

Study of biosynthesis and distribution of minor bioactive compounds in olive oils and tissues applying innovative non-targeted LC-ESI/APCI-MS and GC-APCI-MS approaches.

MP-103 / AN ENERGY-RESOLVED TANDEM MASS SPECTROMETRY APPROACH TOWARDS IDENTIFICATION OF SMALL ISOMERIC METABOLITES

Takemichi Nakamura⁽¹⁾ - Takenori Hama⁽¹⁾ - Masayo Sekimoto⁽¹⁾ - Futoshi Sato⁽²⁾

RIKEN CSRS, Molecular Structure Characterization, Wako, Saitama (1) - Nihon Waters K. K., Market Development Technology Center, Shinagawa, Tokyo (2)

Keywords: Isomer Differentiation, Energy-resolved (ER) MS/MS, CID, Ion Mobility, Collision Cross Section

Introduction

Identification of isomeric organic molecules would become a big challenge for mass spectrometry with ambient in situ setup where no chromatography is available. We explore a viable generic strategy for differentiation of small isomeric metabolites in such situations; the potentials and limitations of widely available low-energy CID MS/MS and ion mobility analysis have been evaluated by using a set of isomeric C4 amino acids as a model.

Methods

Analyte ions were generated by electrospray ionization on directly injected standard solutions. Quadrupole/time-of-flight (Q-TOF) hybrid and triple quadrupole tandem mass spectrometers were used for energy-resolved (ER) CID MS/MS experiments. A traveling wave ion mobility analyzer built in one of the Q-TOF systems was used for ion mobility analysis in N₂. For the CID experiment, He, Ne, Ar, or N₂ was used as the collision gas.

Results

Six C4 amino acids, *i.e.* 2-aminobutyric acid (2-ABA), 3-aminobutyric acid, 4-aminobutyric acid, 2-aminoisobutyric acid (2-AIBA), 3-aminoisobutyric acid, and *N,N*-dimethylglycine (DMG) are one of the simplest examples of small isomers. Among them, 2-ABA, 2-AIBA, and DMG gave similar CID spectra that were attributable to common major fragmentation pathway, loss of 46 Da (CO and H₂O) from protonated molecules (m/z 104) to form m/z 58. Obviously the spectral patterns were condition dependent; however, the difference in energy requirement to give m/z 58 from each isomer could be clearly visualized in breakdown diagrams obtained by ER MS/MS experiments. On the other hand, ion mobility analysis was shown to have limited discrimination power on these isomers. The predicted collision cross section (CCS, calculated by IMoS software [1]) of 2-AIBA protonated molecule was roughly 2% larger than those of the

other isomers whereas the predicted CCS of m/z 58 from 2-ABA was 3% larger than the CCS of m/z 58 from others.

Conclusions

An ER MS/MS approach by using commonly available instruments was shown to be useful for differentiation of small isomeric organic molecules. As the lab-frame collision energy values depend on geometry of the CID device, collision gas species, and gas pressure, method for standardize collision energy axis should be established to facilitate platform independent differentiation and identification of isomers.

Novel Aspect

Potentials and limitations of an ER-MS/MS-based integrated approach for differentiation and identification of small isomeric organic molecules have been evaluated.

References

1. Larriba-AndaluzC., Hogan, C. J. Jr., J. Chem. Phys. 141, 194107 (2014).

MP-104 / HIGHLY EFFICIENT ON-LINE PROTEIN DIGESTION AND PEPTIDE ANALYSIS: EXPLOITING NEW SHORT-BED TRYPSIN-REACTORS

Marina Naldi⁽¹⁾ - **Urh Černigojč**⁽²⁾ - **Ales Štrancar**⁽²⁾ - **Manuela Bartolini**⁽¹⁾

Alma Mater Studiorum University of Bologna, Department of Pharmacy and Biotechnology, Bologna (1) - Bia Separations d. o. o, Ajdovščina (2)

Keywords: Trypsin; Immobilized enzyme reactor; Monoliths; Liquid chromatography-mass spectrometry; Proteomics

Introduction

Trypsin is the gold standard enzyme used in shotgun proteomics for the conversion of proteins into MS-friendly peptides[1]. However, long incubation times and possible covalent modifications induced by digestion conditions are common drawbacks of classic in solution protocols. To achieve high digestion efficiency, automation and a significant decrease in analysis time we considered the development of trypsin-based immobilized enzyme reactors (tryp-IMERs).

Methods

Trypsin-based monolithic tryp-IMERs were obtained by covalent immobilization on short-bedCIMac™ columns(5 mm x 5.2 mm I. D.) exploiting different chemistries and immobilization protocols. IMERs were inserted in an integrated platform featuring analytical columns for peptide analysis and an ESI-Q-TOF mass spectrometer. Optimized procedure was validated in terms of intra- and interday and IMER to IMER reproducibility, ruggedness, stability and sensitivity.

Results

The IMER-based analytical platform allowed the on-line digestion and identification of protein substrates, the separation of peptides and acquisition of MS/MS spectra in an automated manner. In the optimized working conditions no memory effects, nor significant nonspecific interactions were detected. Very short digestion time (90 s), good intra-day and inter-day analysis reproducibility (%RSD in terms of %SQ= \leq 1.2) and good sensitivity down to picomol level were achieved. The high digestion efficiency was confirmed by the analysis of high-molecular weight (HMW) proteins such as a pool of human immunoglobulins G and the HMW fraction of plasma proteins, as challenging real sample. Indeed, efficient protein digestion could be achieved albeit sample complexity, allowing the identification of 281 plasma proteins. Furthermore, the excellent sequence coverages achieved with best performing tryp-IMERs were generally higher than those

obtained by in solution digestion (overnight) or with the commercially available tryp-IMER Poroszyme™.

Conclusions

CIMac columns were suitable supports for the development of highly efficient IMERs. Tryp-IMER performance in terms of digestion of single proteins, protein mixtures and real plasma samples were comparable with that of some capillary and nano-IMERs, without requiring hyphenation to micro- or nano-analytical instrumentation. Automation and high digestion efficiency make the short-bed tryp-IMER appealing new tools for shotgun proteomics and protein analysis.

Novel Aspect

Short-bed CIMac columns were for the first time exploited for the preparation of tryp-IMERs which can be inserted in classic LC-MS systems enhancing automation and throughput.

References

1. Tsiatsiani L., Heck A. J., FEBS J. 282, 2612 (2015).

MP-105 / A COMPARISON OF THE RESULTS OBTAINED FOR THE ANALYSIS OF SOUTH AFRICAN PINOTAGE WINES USING GCxGC-TOFMS WITH HYDROGEN AND HELIUM AS CARRIER GASES

Yvette Naudé⁽¹⁾ - ***Peter Charles Gorst-Allman***⁽²⁾

University of Pretoria, Chemistry, Pretoria (1) - Separation Science Consultant, Separation Science, Pretoria (2)

Keywords: Pinotage wine, carrier gas comparison, GCxGC-TOFMS, principal component analysis, taste profile

Introduction

The relatively low price and availability of helium have made it the default carrier gas for the majority of GC and GC-MS users over the last 40 years. Since the helium shortage of 2012 a number of GC methods have been revalidated for use with hydrogen. This has been driven not only by availability, but also by price, as the cost of helium has doubled in the last 24 months. We compare the analysis of pinotage, a uniquely South African wine, using hydrogen and helium as carrier gas.

Methods

We have analysed seven samples of South African pinotage wines using comprehensive gas chromatography – time of flight mass spectrometry (GCxGC-TOFMS) and both hydrogen and helium as carrier gases. The results were then subjected to principal component analysis (PCA) to investigate the comparability of the components identified using the different carrier gases.

Results

Within the carrier gas groups (*i.e.* helium or hydrogen) results were comparable, but some differences were observed between carrier gas groups. These differences relate primarily to detectability and signal to noise variations for many components, and an attempt has been made to clarify and explain the discrepancies. There are no obvious differences in the mass spectra obtained for the same compounds detected with different carriers, and the mass spectra using both carriers give good quality library matches with spectra contained in commercial mass spectral libraries. By and large the same compounds were found in the same samples analysed with different carriers. In addition the wines were compared by a tasting panel, and an attempt has been made to correlate the wines scoring most highly with their chemical composition and price.

Conclusions

It is likely that helium will be required for GC applications for a long time to come, but as its price continues to rise and it becomes harder to source, laboratories will face continued pressure to pursue alternatives. Hydrogen is the most suitable alternative to helium for many analyses and it is likely to be adopted by GC labs because of its cost effectiveness, availability, and the existence of programmes permitting easy adaptation of existing helium based methods.

Novel Aspect

This is the first time that the compounds providing the organoleptic properties of pinotage wine have been compared using different carrier gases, and GCxGC-TOFMS. The study provides information on analysis of wine using hydrogen as carrier gas, and also highlights several caveats when using this approach.

MP-106 / ANALYSIS OF AN ANTIPARASITIC FRACTION AGAINST GASTROINTESTINAL NEMATODES OF SHEEP FROM LEAVES OF MAYTENUS ILICIFOLIA BY LC-ESI-IT-MSⁿ

Micaela Escobal⁽¹⁾ - *Guillermo de Souza*⁽¹⁾ - *América Mederos*⁽²⁾ - *Cristina Olivaro*⁽¹⁾

Universidad de la República, Espacio de Ciencia y Tecnología Química, Tacuarembó (1) - Instituto Nacional de Investigación Agropecuaria, Programa Nacional de Carne y Lana, Tacuarembó (2)

Keywords: Maytenus ilicifolia Leaves extract, antiparasitic fraction, gastrointestinal nematodes of sheep, LC-ESI-IT-MSⁿ.

Introduction

Maytenus ilicifolia is an important folk medicinal plant, largely found in Paraguay, Uruguay, Argentina, and Brazil. Contributions to human health have been reported [1-3]. A purified fraction from the ethanol extract of Maytenus ilicifolia leaves exhibits antiparasitic activity against gastrointestinal nematodes of sheep. The aim of the present study was to perform a preliminary structural characterization of this bioactive fraction by LC-ESI-IT-MSⁿ.

Methods

The ethanolic extract of M. ilicifolia was subject to a bioguided fractionation on a VLC column and subsequent separation on C18 SPE. The antiparasitic activity was determined against parasites in naturally infected sheep in Uruguay, through egg hatching and inhibition of larval migration in vitro assays. The bioactive fraction was investigated by liquid chromatography/electrospray ionization ion trap multiple-stage mass spectrometry in negative ion mode.

Results

The bioactive fraction was adequately separated by chromatography on a C18 reverse phase UHPLC column using a mobile phase gradient of acetonitrile in water at acidic conditions. The mass spectrometer was set to monitor, select and fragment the most abundant ions [M-H]⁻ in the range m/z 200–1500.

According to their [M-H]⁻ ions, characteristic product ions and retention times, several constituents were identified. Catechin and epicatechin were found. Flavonoid glycosides ranging from 1 to 4 monosaccharide units, having mainly quercetin and kaempferol as aglycones moieties were identified. Tannins were also found, the main ones being composed of catechin/epicatechin and afzelechin/epiafzelechin.

The results of the egg hatching and the larval development assays will be presented at the Conference.

Conclusions

The LC-ESI-IT-MSⁿ is a method quick and easy that can be extended or transferred to other purified bioactive fractions from other plant extracts, allowing a preliminary characterization of its main compounds. In this way the structure of the compounds could be related to the biological activity.

Novel Aspect

This work shows that a purified fraction of ethanolic extract of *Maytenus ilicifolia* rich in flavonoids and condensed tannins present a good antiparasitic activity against gastrointestinal nematode of sheep.

References

1. Souza-Formigoni M. L., Oliviera M. G. M., Monteiro M. G., Silveira-Filho N. G., Braz S., Carlini E. A., Journal of Ethnopharmacology, 34, 21 (1991).
2. Cipriani T. R., Mellinger C. G., Souza L. M., Baggio C. S. F., Marques M. C. A., Gorin P. A. J., Saaki G. L., Iacomini M., Journal of Natural Products, 69, 1018 (2006).
3. Baggio C. H., Freitas C. S., Otofujii G. M., Ciprinani T. R., Souza L. M., Sasaki G. L., Iacomini M., Rieck L., Mesia-Vela S., Marques M. C. A., Journal of Ethnopharmacology, 113, 433 (2007).

MP-107 / FIRST REPORT ABOUT THE STRUCTURE OF AN IMMUNOADJUVANT SAPONIN FRACTION FROM QUILLAJA BRASILIENSIS LEAVES BY LC-ESI-IT-MSⁿ

Federico Wallace⁽¹⁾ - **Zhora Bennadji**⁽²⁾ - **Fernando Ferreira**⁽³⁾ - **Cristina Olivaro**⁽¹⁾

Universidad de la República, Espacio de Ciencia y Tecnología Química, TACUAREMBÓ (1) - Instituto Nacional de Investigación Agropecuaria, Programa Nacional de Investigación en Producción Forestal, TACUAREMBÓ (2) - Universidad de la República, Laboratorio de Carbohidratos y Glicoconjugados, Montevideo (3)

Keywords: Triterpenic saponins, Quillaja Leaves extract, Immunoadjuvant saponin fraction, ESI-MSⁿ, UHPLC.

Introduction

Quillaja brasiliensis is an endemic tree species in Brazil and Uruguay. It has been shown that the aqueous extract and some of the purified saponin fractions obtained from its leaves have immunoadjuvant activity comparable with Quil-A®, the main commercial adjuvant product based on Q. saponaria saponins [1, 2]. The aim of this study was to perform a preliminary structural characterization of the immunoadjuvant fraction of saponins designated as Fraction B (FB).

Methods

The aqueous extract of Q. brasiliensis was subjected to a separation on an SPE column to remove non-saponin components and the fractions containing saponins were reunited. This fraction (FB) was investigated by direct infusion and liquid chromatography/electrospray ionization ion trap multiple-stage mass spectrometry in negative ion mode (DI-ESI-IT-MSⁿ and LC-ESI-IT-MSⁿ). The characterization was performed based on MS² and MS³ experiments of the [M-H]⁻ ions.

Results

The saponins contained in FB were adequately separated by chromatography on a C18 reverse phase UHPLC column using a mobile phase gradient of acetonitrile in water at acidic conditions. The mass spectrometer was set to monitor, select and fragment the most abundant ions [M-H]⁻ in the range *m/z* 1300–2000, using different collision energies. One of them allowed the production of the precursor ion [M-H]⁻ in good yield, while the second energy was set to a higher value in order to produce more product ions that provided more structural information.

According to their [M-H]⁻ ions, characteristic product ions and retention times, twenty seven saponins, were tentatively identified [3]. These saponins feature four different triterpenic aglycones, namely quillaic acid, gypsogenin, phytolaccinic acid and its O-23 acetate. All the saponins are bidesmosidic, bearing oligosaccharides chains linked to the aglycone at C-3 and C-28, the last one containing or not a branched dicarboxylic acid linked or not to a terminal Ara residue.

Conclusions

This is the first detailed report about the structure of *Q. brasiliensis* saponins. Further work will be necessary to fully characterize all the immune stimulating saponins present in the *Q. brasiliensis* species.

Novel Aspect

Q. brasiliensis is an alternative source of immunoadjuvant saponins. The production of Fraction B could scale and become a commercial product comparable with Quil-A®.

References

1. Yendo A. C. A., de Costa F., Cibulski S. P., Teixeira T. F., Colling L. C., Mastrogiovanni M., Soulé S., Roehe P. M., Gosmann G., Ferreira F. A., Fett-Neto A. G., Vaccine, 34, 2305-2311 (2016).
2. Cibulski S. P., Silveira F., Mourglia-Ettlin G., Teixeira T. F., dos Santos H. F., Yendo A. C., de Costa F., Fett-Neto A. G., Gosmann G., Roehe P. M., Comparative Immunology, Microbiology & Infectious Diseases, 45, 1-8 (2016).
3. Wallace F., Bennadji Z., Ferreira F., Olivaro C., Phytochemistry Letters, 20, 228-233 (2017).

MP-108 / A VERSATILE SEMI-PREPARATIVE HPLC PLATFORM, INCLUDING HR MS AND NMR, TO SUPPORT PROCESS DEVELOPMENT IN INDUSTRIAL BIOTECHNOLOGY

Burhan Ozalp

DSM, R&D, Delft, Netherlands

Keywords: Chromatographic fractionation, high resolution mass spectrometry, high field NMR, structure elucidation

Introduction

Identification of the chemical structure of impurities and characterization of their biological activities in the early stage of product development as well as in final product is very important in industrial biotechnology. This is a crucial aspect of “process understanding” and “product specification”. Furthermore, regulatory authorities (ICH, FDA, EMA) have established clear guidelines demanding close monitoring, in-depth investigation and unambiguous confirmation of the chemical structure of impurities in drug products.

Methods

Various approaches can be considered to identify compounds but most of them require relatively pure compounds. While MS has a very high sensitivity, NMR identification in the past required more material. Nowadays NMR instruments have higher sensitivity and approximately 50 µg of pure compound is sufficient for the identification of unknowns, meaning that the composition (identity and relative quantities of elements) and constitution (nature and sequence of bonds), can be determined. Such small amounts can be obtained easily from a semi preparative HPLC platform using 10 to 20 mm i.d. columns.

Results

This study presents a workflow to effectively isolate impurities by using semi preparative HPLC followed by high resolution mass spectrometry and high field NMR analyses for structure elucidation. Furthermore, the workflow can be applied as untargeted approach to characterize complex matrices where the whole chromatogram is fractionated and each fraction separately tested on *e.g.* biological activity and further characterized by HR-MS and NMR. The workflow includes checks on purity and stability of collected fractions and concentration and desalting steps before HR-MS and NMR analysis.

Conclusions

A practical workflow has been developed for the identification of the chemical structure of impurities. The workflow is based on isolation of impurities using a semi preparative HPLC and followed by high resolution

mass spectrometry and high field NMR analyses for structure elucidation. The workflow has been applied successfully for structure elucidation of a number of unknown impurities in various biotech products. Furthermore, it has been also demonstrated that this workflow can be used to characterize samples consisting of complex mixtures.

Novel Aspect

Fast and robust protocol for semi-prep chromatography, HR-MS and NMR workflow for targeted and untargeted fractionation in combination with intermediate checks on purity of collected fractions

MP-109 / VARIOUS CHROMATOGRAPHY APPROACHES HYPHENATED WITH MASS SPECTROMETRY FOR THE DETERMINATION OF FLAVONOIDS AND THEIR METABOLITES

Veronika Pilařová⁽¹⁾ - *Jakub Pavlík*⁽¹⁾ - *Tereza Lacmanová*⁽¹⁾ - *Enrique Espinola*⁽²⁾ - *Lucie Nováková*⁽¹⁾

Faculty of Pharmacy in Hradec Králové, Charles University, Department of Analytical Chemistry, Hradec Králové (1) - University of Granada, University of Granada, Granada (2)

Keywords: ultra-high performance liquid chromatography, hydrophilic interaction chromatography, supercritical fluid chromatography, mass spectrometry, quercetin, phenolic compounds

Introduction

Quercetin is one of the most discussed flavonoid compounds with beneficial effects to human health due to its antioxidant properties. In human body, it is metabolized to two groups of health affecting metabolites with different physico-chemical properties including acidity, molecular weight, and log P. Therefore, it is difficult to determine these compounds in a single analytical run [1][2].

Methods

Three chromatography methods with different selectivity using mass spectrometry (MS) detection in both electrospray positive and negative ion modes were chosen and optimized. These methods included reversed phase ultra-high performance liquid chromatography (RP-UHPLC), hydrophilic interaction chromatography (HILIC), and complementary supercritical fluid chromatography (UHPSFC).

Results

We carried out columns screening with emphasis on separation of isomeric pair of tamarixetin and isorhamnetin. Shield RP18 column using methanol and 0.1% formic acid as the mobile phase enabled the baseline separation of all analytes in RP-UHPLC. Pentahydroxyl, amide, silica, and zwitterion columns with acetonitrile and different aqueous solutions in wide pH range were used for separation in HILIC mode. However, the separation of flavonoids was less successful due to the peak broadening. Ten different columns were screened for SFC. The most promising results were obtained using the diol column and methanol with the addition of organic acids, buffers, ammonium hydroxide, and water. Finally, the MS parameters were tuned to get sensitive detection.

Conclusions

The optimized methods will be compared in the terms of selectivity, retention and sensitivity. The targeted analysis of quercetin and related compounds including small phenolic acids and large metabolites with preserved flavonoid structure in combination with high throughput sample preparation step will be used to explain their metabolism and the health effects.

Novel Aspect

The simultaneous analysis of quercetin and its metabolites with different physico-chemical properties using various separation mechanisms with MS.

References

1. Wang W, Sun C., Mao L., Ma P., Liu F., Yang J., Gao Y., Trends in Food Science and Technology, 56, 21 (2016).
2. D'Andrea G., Fitoterapia, 106, 256 (2015).

The work was financially supported by the project GACR 301/17/05409S and the project STARSS reg. no. : CZ.02.1.01/0.0/0.0/15_003/0000465 funded by EFRR.

MP-110 / MULTIDIMENSIONAL GAS CHROMATOGRAPHY COUPLED TO ISOTOPE RATIO MASS SPECTROMETRY IN SPORTS DRUG TESTING

Thomas Piper - Marlen Putz - Mario Thevis

German Sport University Cologne, Institute of Biochemistry - Center for Preventive Doping Research, Cologne

Keywords: testosterone misuse, carbon isotope ratios, Deans switch, long term detectability

Introduction

Testosterone (T) is frequently misused in sports and its administration is forbidden under the rules of the World Anti-Doping Agency [1]. Its detection remains challenging as T is an endogenous substance. Only by carbon isotope ratio determinations it is possible to distinguish endo- from exogenous T. Sample preparation is commonly based on twofold high performance liquid chromatography clean-up, which is laborious and requires up to 2 working days [2].

Methods

Multidimensional gas chromatography (MDGC) is suitable for the separation of complex matrices. Therefore, one GC equipped with a column of lower polarity is used for sample clean-up and hyphenated to a second GC that employs a more polar column for analyte separation by a pressure-controlled Deans switch device. The effluent from the second GC is partitioned between an isotope ratio and a single quadrupole mass spectrometer to enable simultaneous determinations of isotope ratios and peak purity.

Results

The developed method was fully validated following the guidelines of the World Anti-Doping Agency [3] and employing linear mixing models and reference population derived thresholds [2]. The method demonstrated comparable results to the current screening procedure and sample preparation was shortened to 1 day. Additionally, it was possible to implement measurements of a recently detected new long-term metabolite of testosterone administrations. This epiandrosterone prolongs the detection time by a factor of 2 to 3 [4]. Only for testosterone itself a bias between the screening and the new method was found, which can be partly compensated by population-derived threshold levels.

Conclusions

The developed method employing multi-dimensional gas chromatography instead of high performance liquid chromatography was found to be fit-for-purpose and has been applied to sports drug testing complementing the current routine procedure. Epiandrosterone was found to be a suitable marker for steroid administration and was successfully implemented into doping controls.

Novel Aspect

Multi-dimensional gas chromatography hyphenated to isotope ratio mass spectrometry was successfully employed for doping control analysis substantially reducing the workload for sample preparation.

References

1. World Anti-Doping Agency – Prohibited List 2018. www.wada-ama.org/sites/default/files/prohibited_list_2018_en.pdf
2. Piper T, Thevis M. *Methods in Enzymology*, 596, 403-432 (2017)
3. World Anti-Doping Agency – TD2016IRMS. wada-ama.org/sites/default/files/resources/files/wada-td2016irms-detection_synthetic_forms_eaas_by_irms-en.pdf
4. Piper T, Putz M, Schänzer W, Pop V, McLoed MD, Uduwela DR, Stevenson BL, Thevis M. *Drug Test. Analysis*, 9, 1695-1703 (2017)

MP-111 / MEMS TIME-OF-FLIGHT MASS SPECTROMETER FOR IN SITU GAS MIXTURE ANALYSIS

Frederic Progent⁽¹⁾ - ***Alexandre Sonnette***⁽¹⁾ - ***Jérôme Tupinier***⁽¹⁾ - ***Pierre-Etienne Buthier***⁽¹⁾ - ***Jean-Christophe Lictévout***⁽¹⁾ - ***Sébastien Vigne***⁽¹⁾ - ***Thomas Alava***⁽²⁾

CEA, DAM, DIF, Arpajon, F-91297 (1) - CEA, LETI, MINATEC, Grenoble, F-38054 (2)

Keywords: Miniaturization, μ -TOF, in situ detection, focusing, instrument development

Introduction: Mass Spectrometry detection on-site

Fast and precise in situ gas sensing and identification is of utmost importance in the case of an industrial chemical accident, of atmospheric pollution mapping, for homeland security. The systems used for such kind of purpose need to be of small overall dimensions and weight in addition of their top level performance. Towards this goal, a few research groups [1,2] have developed miniaturized mass spectrometer systems based either on standard or on microfabricated technics.

Methods: Micro Mass spectrometer on chip

This microsystem is fabricated using Micro Electro Mechanical Systems (MEMS) technology. This technic is more promising for size reduction, but also more challenging than standard manufacturing techniques. The MEMS developed in this work have a size of only 1 cm \times 2 cm for a micro-mass spectrometer based on Time-of-Flight (TOF) separation. The micro-TOF ion source is studied and optimized subsystem by subsystem [3,4].

Results: Mass analysis by a linear Micro-TOF

Using electron impact ionization (EI 70 eV), ions are generated in an ionization chamber, then extracted and focalized by electrostatic lenses constituted by 6 electrodes. Ion beam is collected in linear mode in a Micro Channel Plate (MCP) detector and mass spectra are recorded. Using ultrafast pulsing of these ions, the micro ion source is used as a linear TOF spectrometer that is able to separate simple gas mixtures efficiently. Ionization of neutrals atoms is shown to be efficient by electron impact. Using comparison with ion trajectories simulations, the extraction and focalization are optimized to be able to extract the highest ionic current from the ionization chamber. Pure electrostatic simulation can be used as guideline and limits of the simulation values are discussed. The electrode voltage setting of ionization and focusing areas were studied and influences on signal intensity and on mass resolution are followed. Mixes of noble gas, alkane or BTEX were tested in linear micro-TOF.

Conclusions

Miniaturization on chip of mass spectrometer Time-of-Flight is operational and provides repeatable and reproducing results. These devices allow studies of extraction and ion focusing on chip. Electrostatic simulation could be used as global guideline and is a tool for the understanding of ion beam behaviour. Resolution will be improved using orthogonal injection.

Novel Aspect

Gas analyses by linear micro-TOF are available and this work shows encouraging results and pulls the μ -TOF

one step closer towards a fully integrated portable analytical system.

References

1. Syms R., Wright S., J. Micromech. Microeng. 26, 1–28, (2016).
2. Snyder D. T., Pulliam C. J., Ouyang Z., Cooks R. G., Anal. Chem., 88, 2–29, (2016)
3. Vigne S., Alava T., Tasseti C. -M., Duraffourg L., Progent F., Sensors and Actuators B, 243, 690-695 (2017)
4. Vigne S., Alava T., Videlier H., Mahieu R., Tasseti C. -M., Duraffourg L., Progent F., IJMS, 422, 170-176 (2017)

For information please contact: frederic.progent@cea.fr

MP-112 / ATMOSPHERIC PRESSURE CHEMICAL IONIZATION – CONSTRUCTING IONIZATION EFFICIENCY SCALE AND DERIVATIZATION REAGENT SUITABILITY

Riin Rebane - Piia Liigand - Jaanus Liigand - Maarja-Liisa Oldekop - Koit Herodes - Ivo Leito

University of Tartu, Institute of Chemistry, Tartu

Keywords: ionization efficiency; atmospheric pressure chemical ionization (APCI); derivatization reagents.

Introduction

Recent evidence has shown that the atmospheric pressure chemical ionization (APCI) mechanism can be more complex than generally assumed. In order to better understand the factors affecting the APCI, an ionization efficiency scale for an APCI source has been created for the first time. Moreover, based on the scale, five different derivatization reagents for amino acid analysis were tested for APCI suitability (eluent, sensitivity, matrix effects etc).

Methods

The MS responses were recorded in a flow injection mode as previously done by Krueve et al. [1] For every compound, the response versus concentration slope was measured and the ionization efficiencies were expressed in relative terms, using pyrene as the reference compound. Comparison of derivatization reagents was carried out with optimized LC-APCI-MS methods on a reversed phase column. [2]

Results

For the first time an ionization efficiency scale for 40 compounds has been created for APCI source. All logIE values are presented relative to pyrene, chosen as the reference compound with the assigned logIE value 3.41, and range from -0.62 up to 4.32, *i.e.*, spanning almost 5 orders of magnitude. [3]

Suitability of five derivatization reagents for amino acids analysis (diethyl ethoxymethylenemalonate (DEEMM), benzyl ethoxymethylene malonate (DBEMM), dansyl chloride (DNS), 9-fluorenylmethyl chloroformate (FMOC-Cl) and 2,5-dioxopyrrolidin-1-yl *N*-tri(pyrrolidino)phosphoranylideneamino carbamate (FOSF)) were tested for APCI. Acetonitrile proved to be more suitable for the analysis and sample dilution approach in green tea showed that APCI is, in general, less prone to matrix effect than electrospray ionization. FMOC-Cl had poor sensitivity and also high matrix effects (ME). DEEMM and DNS-derivatives were less influenced by the ME than DBEMM and FMOC- derivatives. [2]

Conclusions

The results show that the best ionized compounds in the APCI are not small volatile molecules. Overall the ionization efficiency vs molecular structure trends in APCI are very similar to ESI. results imply that in APCI several ionization mechanisms operate in parallel and mechanism(s) not relying on evaporation of neutral molecules from droplets has significantly higher influence than commonly assumed. Moreover, derivatization reagents for APCI were DEEMM and DNS.

Novel Aspect

For the first time an ionization efficiency scale spanning over 5 logIE units and including 40 diverse compounds has been created for APCI source.

References

1. Krueve A., Kaupmees K., Liigand J., Leito I. Analytical Chemistry, 86, 4822(2014).
2. Oldekop M. -L., Herodes K., Rebane R., International Journal of Mass Spectrometry, 421, 189 (2017).
3. Rebane R., Krueve A., Liigand P., Liigand J., Herodes K., Leito I. Analytical Chemistry, 88, 3435 (2016).

MP-113 / MICROFLUIDIC CE-MS FOR NEWBORN SCREENING: A SINGLE SYSTEM FOR MONITORING SMALL MOLECULE AND PROTEIN BIOMARKERS

Erin Redman - J. Scott Mellors

908 Devices Inc., Life Sciences, Boston

Introduction

After birth several assays are typically performed to determine if an infant suffers from conditions that could have long-term health effects. The compounds measured range from small molecules to larger proteins. Few techniques are capable of effectively analyzing both small and large molecules rapidly and simply. In this work we present a microfluidic CE-MS system for analyzing both small molecules and intact proteins related to newborn screening assays.

Methods

Microfluidic chips with a 10 cm separation channel were used for small molecule assays and chips with a 22 cm separation channel were used to analyze protein biomarkers (908 Devices Inc.). Samples for small molecule assays a liquid extraction was performed to elute the compounds from dried blood spots. Protein biomarkers were measured from whole blood or isolated erythrocytes. MS analysis was performed using a Thermo Orbitrap MS.

Results

Results will be presented demonstrating microfluidic CE-MS analysis of several common biomarkers regularly monitored during new born screening. Using the same method, amino acids and acyl carnitines can be monitored to assess whether samples are indicative of Phenylketonuria (PKU) or fatty acid metabolic disorders, such as Carnitine acylcarnitine translocase deficiency (CACT). The small molecule biomarkers are able to be quantified over the appropriate physiological ranges indicative of disease states. Hemoglobin subunits from whole blood or purified erythrocyte lysates were analyzed to detect hemoglobin variants indicative of hemoglobinopathies, such as Sickle Cell Anemia. The alpha, beta, gamma, and abnormal subunits were resolved in the microfluidic CE separation within 4 minutes and the identity of the subunits was confirmed by mass. In the samples analyzed, beta chain variants indicative of sickle-cell disease were characterized.

Conclusions

The microfluidic CE-MS platform can be used to perform rapid small molecule and protein assays related to new born screening. Additionally, multiple small molecules can be assessed with a single method providing the ability to screen for different disease states in a single run. By analyzing hemoglobin at the subunit level, minor modifications can be identified and detected by both their mass and anticipated change in migration times due to modification.

Novel Aspect

In this work a single system is used to demonstrate new born screening assays for small molecule and protein biomarkers.

MP-114 / DEVELOPMENT OF RAPID SCREENING APPARATUS WITH DIP AND PI/MS FOR ADDITIVES IN POLYMER MATERIALS

Yuji Mishima ⁽¹⁾ - Takahisa Tsugoshi ⁽²⁾

KMTL, KMTL, Tokyo, Japan (1) - NMIJ, AIST, Tsukuba, Japan (2)

Keywords: direct inlet probe, fragmentless ionization, screening, silent change, additive,

Introduction

The globalization of the economy has made risks of quality control due to the “Silent Change” apparent in manufacturers, particularly in the electrical and electronic industries. The Silent Change refers to material specifications being modified without approval in the supply chain, and leads sometimes to serious accidents, that causes to product recall and the loss of the brand image.

For its quality control, rapid screening technique is required with all detection of additive species because of many specimens of lots, and with simple operation of apparatus.

Methods

For all detection of additives, it is better to use DIP, direct inlet probe that extract additives thermally in vacuum. However low boiling point chemicals are vaporized and disappeared at measurement, so that technique to cover the specimen is used, named SET, Suppressing Evaporation Technique. Separation technique as GC and LC takes time, so that it is required a fragmentless ionization that realizes informational separation in a mass spectrum because one peak indicates one species.

Results

JMS-Q1500GC by JEOL was accepted as a base instrument for the developed new apparatus. This instrument has both options of DIP and photoionization that is one of the fragmentless ionization techniques. However these options are exclusive. To be consistent with both, an alteration had been required. Further new design of DIP had been required. We manufactured a new design vacuum chamber and DIP by way of trial. At this stage of development, performance of the trial apparatus had been checked with some specimens. For example, good calibration curve has been obtained from Vitamin E within 10 min., but not from Decabromodiphenyl because of the design of instrument.

In this design, DIP was inserted to an underside of the chamber, so that the operability was not good. The performance of coupling of DIP and Photoionization had already ascertained, consequently we moved to the next step to solve the operability problem. Newer designs of DIP and vacuum chamber were required for that DIP is possible to insert to a front of the chamber. And the new design has resolved the detection problem of Decabromodiphenyl.

Conclusions

The new apparatus with both of DIP and photoionization consistently has been successfully developed. The measurement is without property dependence from chemicals, *e.g.*, a difference of boiling point, as all detection.

Novel Aspect

The developed apparatus offers high level quality control to avoid the silent change with rapidity and easy operation. It can be a deterrent to the silent change.

MP-115 / COMPARISON OF THE DISPLACEMENT MODE VERSUS GRADIENT MODE FOR PROTEOMICS USING AN ONLINE 2D-LC-MSMS SYSTEM

Hartmut Schlüter

University Medical Center Hamburg-Eppendorf, Mass Spectrometric Proteomics, Hamburg

Keywords: on-line-two-dimensional-liquid-chromatography-tandem-mass-spectrometry (2D-LC-MSMS); displacement chromatography; proteomics;

Introduction

The extreme complexity of proteomes is a major challenge for LC-MSMS-based proteomics. Two-dimensional-(2D)-LC applying cation-exchange chromatography (CEX) in the first dimension has been successfully used for increasing the number of identified proteins. In this study it was investigated if CEX performed in the displacement mode (DM) within an online-2D-LC-MSMS system offers advantages in comparison to the gradient mode (GM) for the analysis of proteomes.

Methods

Tryptic peptides of Hela cells (5 µg l total) were applied to an online-2D-LC-MSMS system (1st dimension: CEX column; 2nd dimension reversed-phase column). Tryptic peptides were eluted from the CEX either by pulsed injections of increasing amounts of ammonium acetate (GM) or by pulsed injections of spermine (DM). Resulting MS data were subjected to OpenMS, MaxQuant and mathematica scripts.

Results

The total number of identified peptides and proteins as well as the sequence coverage was significantly higher by using DM compared to GM. A major reason for this phenomenon is, that DM was giving a significantly better separation via CEX of doubly charged tryptic peptides (charge state in the mobile phase of the LC), which represent the majority of tryptic peptides in mammalian proteomes. The results are in good agreement with former studies about the application of the displacement mode for the separation of complex mixtures in combination with mass spectrometry [1, 2]. However, in comparison to the former study [2] in the current study a big step forward was achieved regarding the operating time (reduction from 6 days to approximately 12 hours) of the 2D-LC-MSMS analysis. In addition the total sample amount was reduced down to 5 µg.

Conclusions

In comparison to gradient mode the displacement mode for cation-exchange chromatography within an online system comprising two-dimensional liquid-chromatography coupled to tandem mass spectrometry for identifying proteomes via tryptic peptides is significantly more efficient.

Novel Aspect

First study comparing the efficacy of displacement mode versus gradient mode for proteomics using online-2D-LC-MSMS

References

1. Ahrends R., Lichtner B., Bertsch A., Kohlbacher O., Hildebrand D, Trusch M, Schlüter H. J Chromatogr A. 1217:3321-9 (2010).
2. Trusch M., Tillack K., Kwiatkowski M., Bertsch A., Ahrends R., Kohlbacher O., Martin R., Sospedra M., Schlüter H., J Chromatogr A. 1232, 288-94 (2012).

For information please contact: hschluet@uke.de

MP-116 / LOW-FLOW ATMOSPHERIC PRESSURE CHEMICAL IONIZATION BASED ON A GAS DYNAMIC VIRTUAL NOZZLE NEBULIZER

Timotej Strmeň - Vladimír Vrkoslav - Ondřej Pačes - Josef Cvačka

Institute of Organic Chemistry and Biochemistry of the Czech Academy of Sciences, Department of Analytical Chemistry, Faculty of Science, Charles University in Prague, Prague

Keywords: low-flow APCI, GDVN

Introduction

Atmospheric pressure chemical ionization (APCI) is a useful tool for analyzing small and relatively stable molecules, *e.g.* lipids [1], at milliliter-per-minute sample flow rates. Despite the previous effort to develop low-flow APCI [2,3], no such an ion source has been commercialized yet. We are currently developing a nebulizer for an APCI source operable at sub-microliters per minute flow rates.

Methods

The sprayer was manufactured from a borosilicate glass tube (I.D. 1.15 mm) and two fused-silica capillaries, which were fixed inside the tube: (i) tapered-tip capillary for the sample and (ii) capillary for gas delivery. The device generated liquid jet thanks to gas-dynamic forces exerted by a co-flowing nitrogen (gas dynamic virtual nozzle, GDVN [4]). It served as a nebulizer for low-flow APCI.

Results

The spray created by the GDVN nebulizer operated at 120 ml/min nebulizing gas flow rate and 200 nl/min liquid sample flow rate had a conical shape with the diameter of 1.5 mm at the distance of 2.4 cm. To facilitate fast solvent evaporation, a short, resistively heated glass tube was placed between the sprayer orifice and corona discharge region. Test compounds including acridine in methanol and acetonitrile, cholesterol in toluene, cis-11-hexadecenol in acetonitrile (all in positive mode), and 5-nitrosalicylic acid in methanol (in negative mode) were successfully detected. Their spectra were similar to those recorded from the conventional APCI. The signal stability was rather moderate (RSDs of acridine signal recorded for ten minutes were about 30 – 40 %). High signal intensities of the analytes suggested the enhanced sensitivity of the GDVN APCI source.

Conclusions

An APCI source with a GDVN nebulizer was found appropriate for nanoliter-per-minute sample flow rates. Various compounds were detected using a liquid flow rates down to 200 nl/min. No clogging of the sprayer was observed throughout the measurements, which, together with its high sensitivity, makes it a good candidate for a nano-HPLC and mass spectrometry hyphenation.

Novel Aspect

For the first time, GDVN sprayer was successfully used as a nebulizer for liquid samples in an ion source.

References

1. Byrdwell W. C. :Lipids, 36, p. 327 – 346 (2001)
2. Takeda S., Tanaka Y., Yamane M., Siroma Z., Wakida S., Otsuka K., Terabe S. :J Chromatogr A, 924, p. 415-420 (2001)
3. Östman P., Marttila S. J., Kotiaho T., Franssila S., Kostianen R. :Anal Chem, 76, p. 6659-6664 (2004)
4. DePonte D. P., Weierstall U., Schmidt K., Warner J., Starodub D., Spence J. C. H., Doak R. B. :J Phys D:Appl Phys, 41, 195505 (2008)

MP-117 / ANALYSIS OF ENZYMES IN LAUNDRY DETERGENT USING SDS-PAGE AND MALDI-TOF-MS

Zahra Talebpour⁽¹⁾ - **Ghazaleh Abedi**⁽¹⁾ - **Alireza Ghassempour**⁽²⁾

Department of Chemistry, Alzahra University, Tehran (1) - Faculty of Phytochemistry, Shahid Beheshti University, Tehran (2)

Keywords: Enzyme, laundry detergent, SDS-PAGE, MALDI-TOF

Introduction

There is a major challenge for the development of method in order to identify the enzymes in raw and manufactured materials in cosmetic industries[1, 2]. In this study we used SDS-PAGE and MALDI-TOF-MS instrument for identification of main types of enzymes used in laundry detergents: protease, lipase, cellulose, amylase, and mannanase.

Methods

A certain weight of granular solid concentrates, was dissolved in Tris-base, and the suspension was stirred for 30 min. After centrifugation, equal volume of TCA and supernatant were mixed, and enzyme precipitation was accomplished with TCA protocol. These mixtures were separated according to their molecular weights on 12% Bis-Tris gels and visualisation was performed using Coomassie brilliant blue. After cutting of gel and remove of Coomassie reagent, the mass values has obtained by using MALDI-TOF-MS.

Results

Samples were load on the gel, and the separation was done based on the peptides' molecular weight. The main protein bonds were detected at approximately 26, 72& 40, 32, 55, 55, and 35 kDa for Everlase, Carezyme, Mannaway, Termamyl, Duramyle, and Lipolase respectively. In order to obtain precise mass values, the same solutions were analyzed by MALDI-TOF-MS. Finally, the precise mass values were obtained, 26.6, 72.5& 35.6, 33.7, 55.6, 55.6, and 31.7 kDa for Everlase, Carezyme, Mannaway, Termamyl, Duramyle, and Lipolase respectively and results have compared with reported MS results of these enzymes. Also, the protein bonds were detected successfully in ternary, quaternary and quinary mixtures of crude enzymes. As a final point the granular enzymes were analyzed with SDS-PAGE and MALDI-TOF. According to the obtained results, we can identify every kind of granular enzymes.

Conclusions

In this study, determination of protease, lipase, cellulose, amylase, and mannanase was accomplished by SDS-PAGE and MALDI-TOF-MS. The MALDI-TOF-MS as a high-throughput screening method can be proposed for monitoring of enzymes used in laundry detergents, owing to its advantages such as speed, sensitivity, reproducibility and label-free readout.

Novel Aspect

SDS-PAGE and MALDI-TOF-MS can be proposed as a versatile quality control approaches for detergent industry to identify used enzymes.

References

1. Gaubert, A., Jeudy, J., Rougemont, B., Bordes, C., Lemoine, J., Casabianca, H., & Salvador, A., Analytical and bioanalytical chemistry, 408:17, 4669-4681, (2016).
2. Beneito-Cambra, M., Bernabé-Zafón, V., Herrero-Martínez, J. M., Simó-Alfonso, E. F. and Ramis-Ramos, G., Talanta, 79:2, 275-279, (2009).

MP-118 / MECHANISTIC REACTION INSIGHTS BY *IN-SITU* MONITORING OF AN AACVD REACTION USING GC-MS: STUDYING THE ZnFe₂O₄ THIN FILM GROWTH PROCESS

Ime C. Usen - Upul Wijayantha-Kahagala-Gamage - James C. Reynolds

Department of Chemistry, Loughborough University, Loughborough, United Kingdom

Keywords: AACVD, Gas Chromatography Mass Spectrometry, Reaction Monitoring, Volatile profile, Exhaust gas

Introduction

A variant of the conventional chemical vapour deposition (CVD) process, the aerosol-assisted (AA) CVD method, uses aerosol droplets to transport the precursors, with the aid of inert or reactive carrier gases to the substrate for thin film coating. It is important to monitor the reaction that takes place in the AACVD system. A novel adaptation of GC-MS is developed which enables direct monitoring of AACVD exhaust gases for the first time providing valuable information.

Methods

The AACVD setup described by Tahir et al. (1) was directly interfaced via a heated transfer line and semi-automated port valve to GC-MS with DB-1 column. AACVD was conducted by systematically increasing the temperature from 300 to 550°C in 50°C intervals. Measurements were performed under programmed temperature conditions and the resulting films were characterised by X-ray diffraction and Raman spectroscopy to understand material properties and chemical nature.

Results

The exhaust gases were monitored for the solvent (ethanol) and the individual precursor components; zinc acetate and iron (III) chloride and $[\text{Fe}_2(\text{acac})_4(\text{dmaeH})_2][\text{ZnCl}_4]$ the precursor. The GC-MS peaks were successfully observed for the entire experiment at the various temperatures. Key observation noted were: the exhaust gas profile changes depending on temperature; the ethanol changes into acetaldehyde if no precursor is present; when the precursor dissolves in ethanol, other species are observed which are from either the direct breakdown of the precursor, or reaction in the AACVD system, acetone, ethyl ether, ethyl acetate, butenal, 2-butanone and 3-penten-2-one.

The optimal temperature required to produce ZnFe_2O_4 was 450°C from precursor in ethanol. The exhaust gas profile obtain from the GC-MS was confirmed through (XRD) and Raman spectroscopy studies. This study shows that this method can be used to give a mechanistic insight into the processes that takes place in the AACVD system for the first time, with potential for further research.

Conclusions

The volatile profile obtained from the exhaust of an AACVD system provides clear insight to the progression of the overall reaction. Observation of specific reaction products acetone, ethyl acetate, ethyl ether and butenal in exhaust gas phase can be used as an indication for the formation of ZnFe_2O_4 film from precursor $[\text{Fe}_2(\text{acac})_4(\text{dmaeH})_2][\text{ZnCl}_4]$ in ethanol at 450°C. The technique has the potential of being incorporated into simple sensor technology for a wider use.

Novel Aspect

A new application of GC-MS has been developed for in-situ monitoring of AACVD process, with mechanistic information from the volatile profile of the exhaust gas being so vital.

References

1. Tahir A. A., Wijayantha K. G. U., Mazhar M., McKee V., Thin Solid Films. 518(14), 3664. (2010).

MP-119 / PULSED EXTRACTION MODE FOR OFF-AXIS TIME OF FLIGHT MASS SPECTROMETER

Guanyi Wei - Zhiming Li - Lihua Zhai - Xianglong Yuan - Hu Deng - Xiaopan Shen - Wenliang Wang

Northwest Institute of Nuclear Technology, Northwest Institute of Nuclear Technology, Xi'an

Keywords: Sensitivity, Time of flight, Ion Extraction, Pulse

Introduction

Better sensitivity is the long-lasting pursuit of an instrument designer. In most cases of the time of flight (TOF) mass spectrometer (MS), ions were extracted from the ion source continuously and path through the

push area of the TOF-MS. Hereafter, most of the ions in the push area was repeatedly accelerated to the free field area, and will finally reach the detector in sequence, which is decided by their mass to charge ratio.

Methods

A pulsed extracting field and an accelerating field was introduced into a TOF-MS between the ion source and the push area. Thus the ions could be extracted and spread along its path to the push area. Where an ion's position among them depends on its mass to charge ratio. If the push pulse is synchronized to a certain time, it could seize the interested ions by their mass to charge ratio, leave other ions passing by.

Results

Experiments were done on an off-axis TOF with thermal ionization ion source and Cs sample with contamination of Na and K. Working in continuously extracting mode, the acquired mass spectrum will show that the mass peak of Na is much higher than that of Cs. While in pulsed ion extraction mode and synchronized push, the mass peak of Cs is enhanced to several times higher if Cs was focused. The mass peak of Na could be suppressed to below 20% at the same time. If it is synchronized to Na, the mass peak of Na would be enhanced and Cs be suppressed. In this experiment, the difference of the synchronize time, or push pulse delay, between focused to Na and Cs is about 5 micro-seconds, which can be easily handled.

Conclusions

Pulsed ion extracting mode could enhance the interested mass peak by several times and greatly suppress other mass peaks at the same time. Comparing to the continuously extracting mode, it could greatly improve the sensitivity of the instrument without significant change for an off-axis TOF-MS.

Novel Aspect

Pulsed ion extracting mode could enhance the interested mass peak by several times, which shows its great talent on improving the sensitivity of an off-axis TOF-MS, and other instruments.

This work is financially supported by the National Key Scientific Instrument & Equipment Development Projects(Grant No. 2012YQ250003).

MP-120 / PROBING THE LIMITS OF ON-LINE REACTION MONITORING BY SMALL FOOTPRINT ESI-MS AND NMR

Christopher Whitmore⁽¹⁾ - *Alan Kenwright*⁽¹⁾ - *Andrew Ray*⁽²⁾ - *Steven Coombes*⁽²⁾ - *Jackie Mosely*⁽¹⁾

Durham University, Department of Chemistry, Durham (1) - AstraZeneca, Global Chemical Development, Pharmaceutical Development, Macclesfield (2)

Keywords: on-line, reaction monitoring, ESI-MS, compact

Introduction

The proliferation of relatively cheap small-footprint MS instruments presents an opportunity for their application to on-line reaction monitoring, particularly when used in combination with other small footprint analytical instruments [1]. The opportunities and limitations of small footprint ESI-MS in reaction monitoring are explored here.

Methods

A combined system consisting of a Waters QDa single quadrupole ESI-MS instrument and a Magritek Spinsolve 43 MHz NMR instrument were used to monitor a variety of reactions, with on-line dilution provided by a system of Teledyne-SSI LS class HPLC pumps and a Rheodyne Mass Rate Attenuator.

Results

The versatility of the combined ESI-MS/NMR system are demonstrated, as well as the value of a combinatorial approach to reaction monitoring, with examples presented in which data from a single technique would be misleading, and in which the combined data can be used to overcome inherent issues with the techniques in question. The capabilities of the system with regard to extremes of temperature, pH and other reaction conditions are reported.

Conclusions

System is shown to be robust to a variety of conditions, though the performance of rare-earth magnet based NMR to reactions with off-gassing or extremes of temperature is found to be poor. Challenges are found with the use of the single quadrupole MS instrument for the monitoring of reactions with particularly complex mixtures, with methods developed to correct for suppression in ESI-MS.

Novel Aspect

Combined portable small footprint NMR/MS system utilized on a variety of reactions, capabilities of system explored and discussed.

References

1. Ray, A. ; Bristow, T. ; Whitmore, C. ; Mosely, J., Mass Spec. Revs., doi 10.1002/mas. 21539

MP-121 / IMPROVEMENT OF ELECTROSTATIC QUADRUPOLE ZOOM LENS FOR FIXED MULTI-COLLECTOR SYSTEM IN LRI-MS

Xianglong Yuan - Zhiming Li - Lihua Zhai - Jiang Xu - Xiaopan Shen - Hu Deng - Wei Wang

Northwest Institute of Nuclear Technology, Northwest Institute of Nuclear Technology, Xi'an

Keywords: LRIMS; Multi-collector; Zoom lens; Electrostatic quadrupole

Introduction

The Laser-Resonance-Ionization Mass Spectrometer (LRI-MS) is a type of mass spectrometer with high elemental selectivity and high sensitivity, which is able to overcome the possible isobaric interferences [1]. Although frequently utilized in the reported LRI-MS facilities or instruments [2-5], time-of-flight (TOF) analyzer is not very appropriate for precise isotope ratio measurements due to its inherited shortcomings of the detection system, such as variation of signal detection efficiency [6]. In our laboratory, a LRI-MS system equipped with a magnetic analyzer (MA) has been used to measure the isotope ratios of plutonium [7] and tin [8] with peak-jumping method. Recently, a LRI-MS system equipped with an electrostatic analyzer (ESA), a magnetic analyzer (MA) and a fixed multi-collector system has been developed for precise isotope ratio measurement. In this instrument, the distance between the ion beams of different mass numbers at the focal plane is proportional to the relative mass dispersion. In order to achieve simultaneous detection of ions with different central mass by using the fixed multi-collector system, a zoom lens based on electrostatic quadrupoles has been developed to adjust the distances of the ion beams at the focal plane [9]. However, the zoom effect, that is, the horizontal magnification simulated by SIMION is quite different from that predicted by theoretical model.

In this work, the reason why the zoom effect in the simulation is so different from that in theoretical model is carefully analyzed. Based on the analysis, the zoom lens was improved to obtain a more accurate zoom effect. By using a home-made diagnostic tool capable of visualizing the shape of ion beam, the zoom effect of the improved zoom lens was evaluated experimentally.

Theory and Methods

Quadrupole lenses are focusing in one plane and defocusing in another perpendicular one [10]. Thus,

several quadrupole lenses could be combined for a useful lens system. A perfect quadrupole consists of four hyperbolic rods. Since hyperbolic rods are difficult to be manufactured and fabricated, round rods with a diameter of 1.148 times the desired diameter of the quadrupole field were used instead [11]. The shortcoming of such quadrupole lenses is that it occupies more space. In electrostatic field, the solution, which satisfies Laplace equation and given boundary conditions, is unique. Thus, a rectangular boundary field which saves more space could be used. Several segmented electrodes can be used to fit the rectangular boundary of the quadrupole field. However, high-order field components will be introduced. In order to reduce the high-order field components, the structural parameters of the electrodes must be carefully designed. In ideal quadrupole lenses, the electrostatic field strengths are assumed to start and end abruptly, which is impossible in real lenses. To eliminate the influence of the fringing field, an effective length should be evaluated accurately. One method to evaluate the effective length of the quadrupole with round rods is to use the empirical formula $l_e = l_0 + 1.1 a_0$, where, l_e , l_0 and a_0 mean the effective length, real length and aperture of the quadrupole, respectively[12]. But it is not verified in the segmented quadrupole with a rectangular boundary. Another method to evaluate the effective length is fringing field integration. It needs the accurate fringing field distribution, which can be determined by the finite element method (FEM) or finite difference method (FDM). To achieve point-to-point focusing in horizontal direction and make the lens system not too complicated, a better choice is to use two quadrupoles, the parameters of which can be calculated based on the design requirement and ion optics theory.

Experiment and Results

In previous work [9], 24 electrodes were used to fit the rectangular(30 mm×30 mm) field boundary in which 11×2 electrodes were installed in vertical position while 1×2 electrode in horizontal position. It was found that high-order field components and a large constant field existed, the potential of which was 18.5% of the quadrupole field amplitude voltage. To eliminate the constant field and reduce the high-order field components, the rectangular was changed into 60 mm×30 mm, and the left and right electrodes were replaced with three electrodes each side, the potentials of which were equal and the positions of which were hyperbolic. This improvement was verified by the simulation results. After simulated the fringing field distribution by the finite element method, it was observed that the effective length of the field was roughly equal to the actual length of the electrode plus 25 mm.

The improved zoom lens was manufactured and fabricated into the LRIMS. By visualizing the ion beams of europium and barium isotopes at the focal plane, respectively, using a home-made diagnostic instrument, the zoom effect was evaluated experimentally. The results showed that the horizontal magnification ratio is between 0.6 and 1.2 when the longitudinal magnification ratio varies from 0.7 to 1.1. The relative difference between the measured value of the magnification ratio and the predicted value is better than 1%. By fine adjustment of the real voltage of the electrodes, relative to the theoretical voltage, simultaneous detection of the isotopes of different elements has been realized by the fixed multi-collector system.

Conclusions and Novel Aspect

Based on the simultaneous detection requirement of the LRI-MS equipped with an ESA, a MA and a fixed multi-collector system, an ion zoom lens based on the quadrupoles was developed and improved. The influence of the fringing field of the quadrupole is evaluated. The experimental results showed that the lens system achieved accurate zoom effect of the ion beams in the range of 0.6 to 1.2.

References

1. J. D. Fassett, L. J. Moore, J. C. Travis, et al. Laser resonance ionization mass spectrometry [J]. *Science*, 1985,230(4723) :262-267.
2. I. V. Veryovkin, W. F. Calaway, J. F. Moore, et al. A new time-of-flight instrument for quantitative surface analysis [J]. *Nuclear Instruments and methods in Physics Research*, 2004, 219-220(7):473-479. doi:10.1016/j.nimb.2004.01.105
3. C. Gruning, G. Huber, P. Klopp, et al. Resonance ionization mass spectrometry for ultratrace analysis of plutonium with a new solid state laser system [J]. *International Journal of Mass Spectrometry*, 2004,

235:171-178. doi:10.1016/j.ijms.2004.04.013

4. S. A. Crowther, R. K. Mohapatra, G. Turner, et al. Characteristics and applications of RELAX, an ultrasensitive resonance ionization mass spectrometer for xenon [J]. *Journal of Analytical Atomic Spectrometry*, 2008, 7(7):938-947. doi:10.1039/B802899K
5. Yoshihiro Iwata, Chikara Ito, Hideki Harano, et al. Performance evaluation of a resonance ionization mass spectrometer developed for the FFDL system of fast reactors [J]. *Journal of Nuclear Science and Technology*, 2014, 51(4):465-475. doi:10.1080/00223131.2014.874960
6. A. Rowland and J. A. Holcombe. Evaluation and correction of isotope ratio inaccuracy on inductively coupled plasma time-of-flight mass spectrometry [J]. *Spectrochimica Acta Part B Atomic Spectroscopy*, 2009, 64,35-41. doi:10.1016/j.sab.2008.10.030
7. SHEN Xiao-pan, LI Zhi-ming, ZHAI Li-hua, et al. Selective Ionization of Trace Plutonium from Uranium-Plutonium Mixture by Laser Resonance Ionization Mass Spectrometry[J]. *Journal of Chinese Mass Spectrometry Society*, 2013, 34(5): 269-273. doi: 10.7538/zpxb.2013.34.05.0269
8. SHEN Xiao-pan, LI Zhi-ming, Wang Wen-liang, et al. Measurement of Tin Isotope Ratios by Laser Resonance Ionization Mass Spectrometry[J]. *Chinese Journal of Analytical Chemistry*, 2017, 45(3): 342-347. doi: 10.11895/j.issn.0253-3820.160840
9. Yuchi Chen, LihuaZhai, Hu Deng et al. The design of electrostatic quadrupole zoom lens for fixedmulti-collector system at mass spectrometer. *International Journal of Mass Spectrometry*. 2016 (408) 1–8. doi:10.1016/j.ijms.2016.09.001
10. Hermann Wollnik, *Optics of charged particles*, Academic Press, Inc.1987,48
11. Raymond E. March, John F. J. Todd, *Quadrupole ion trap mass spectrometry-2nd ed.*, Pulished by John Wiley & Sons, Inc. 2005, 39
12. F. Grime, G. W. Watt, *Beam Optics of Quadrupole Probe-forming Systems*, Adam Hilger, 1984.

MPS-S05 - Volatile molecules & nose-space analysis

Chairs: Egmont Rohwer, Giorgio Mellerio

MP-122 / OVERVIEW OF THE VOLATILE COMPONENTS OF THE CORSICAN LIVERWORT SCAPANIA UNDULATA

Anaïs Pannequin⁽¹⁾ - **Aura Tintaru**⁽²⁾ - **Jean-Marie Desjobert**⁽¹⁾ - **Jean Costa**⁽¹⁾ - **Alain Muselli**⁽¹⁾

Université de Corse, UMR CNRS 6134, Laboratoire Chimie des Produits Naturels, BP 52, 20250 Corte, France
(1) - Aix Marseille Univ, CNRS, Institut de Chimie Radicalaire, UMR 7273, 13397 Marseille, France (2)

Keywords: Scapania undulata, essential oil, hydrosol, SPME

Introduction

S. undulata is a liverwort living wild on the rocks of deciduous broadleaved forests of Corsica [1]. To our knowledge, four studies deals with the chemical composition of *S. undulata* essential oil from Europe [2–5]. Some sesquiterpenes such as longiborane and himachalane derivatives were identified as main components. The aim of the present work was to determine and characterize the volatile compounds of *S. undulata* from Corsica using hyphenated methods.

Methods

Plant material was collected in several Corsican forests and dried at room temperature. Essential oils and hydrosols were prepared by hydrodistillation as well as the volatile fraction was extracted by HS-SPME (HeadSpace - Solid Phase MicroExtraction) under optimized conditions. Et₂O extract was prepared by

maceration from powdered plant material. The plant volatiles were investigated using an analytical procedure including LC, LLE, GC-FID, GC-MS, ESI-MS, tandem MS and NMR.

Results

The present study focuses on the volatile components of *S. undulata* from Corsica. The chemical investigations of the four extraction matrices allow determining different chemical profiles. Essential oil (EO), hydrosol (HY) and Et₂O extract (EXT) was dominated by oxygenated sesquiterpenes (EO 63.5%; HY 82.0% and EXT 63.9%) while the volatile fraction (VF) extract by HS-SPME was dominated by hydrocarbon compounds (VF 61.2%). The main components of essential oil and hydrosol of *S. undulata* were respectively epicubenol (49.0% and 64.5%), β -Muurolol (5.0% and 7.2%) and longiborneol (3.0% and 7.2%) while epicubenol (44.4%) with longipinanol (11.9%) were the main components of Et₂O extract. As previously reported, longipinanol is rather labile and was almost completely degraded during hydrodistillation, but was detected in considerable amounts in the ether extract[3]. In addition, epicubenol (27%), calarene (9.6%), were the main components of the volatile fraction sampled by HS-SMPE. As epicubenol, himachalol, longiborneol and logipinanol were absents to our spectral library, their identifications were ensured by NMR and ESI-MS experiments.

Conclusions

The study provides an overview of the volatile components of *S. undulata*. Sesquiterpene compounds were the dominant natural products of the plant extracts obtained by hydrodistillation, HS-SPME and Et₂O maceration. Oxygenated compounds characterized essential oil, hydrosol and Et₂O extract while the volatiles emitted by the plant were largely hydrocarbons. The study confirms the necessity to use different sample preparation methods in order to study the complete chemical diversity of the plant.

Novel Aspect

The study reports for the first time the volatile components of *S. undulata* from Corsica. To our knowledge, the report of epicubenol as main essential oil component was original, as well as, the study of the volatile components emitted by *S. undulata*. HS-SPME analysis allows approaching of chemicals emitted by the plant and the identification provides information concerning the allelopathic properties of the liverwort.

References

1. Sotiaux A, Sotiaux O, Vanderpoorten A. Additions to the bryophyte Flora of Corsica. 2008; 29: 267–74.
2. Adio AM, Paul C, Kloth P, König WA. Sesquiterpenes of the liverwort *Scapania undulata*. *Phytochemistry* 2004; 65: 199–206.
3. Andersen NH, Bissonette P, Liu C-B et al. Sesquiterpenes of nine european liverworts from the genera, *Anastrepta*, *bazzania*, *jungermannia*, *lepidozia* and *Scapania*. *Phytochemistry* 1977; 16: 1731–51.
4. Connolly JD, Phillips WR, Huneck S. (+)-Ent-epicubenol from the liverwort *Scapania undulata*. *Phytochemistry* 1982; 21: 233–4.
5. König WA, Rieck A, Hardt I, Gehrcke B, Kubeczka K-H, Muhle H. Enantiomeric composition of the chiral constituents of essential oils. Part 2: Sesquiterpene hydrocarbons. *Journal of High Resolution Chromatography* 1994; 17: 315–20.

MP-123 / EFFECT OF GERMINATION TIME, ROASTING LEVEL AND ALTITUDE IN THE AROMA RELEASED BY COFFEE POWDER MONITORED BY PROTON TRANSFER REACTION MASS SPECTROMETRY AND GAS CHROMATOGRAPHY MASS SPECTROMETRY

Martina Bodner⁽¹⁾ - *Ksenia Morozova*⁽¹⁾ - *Prapaijit Kruathongsri*⁽²⁾ - *Prodpran Thakeow*⁽²⁾ - *Sergio Angeli*⁽¹⁾ - *Matteo Scampicchio*⁽¹⁾

Free University of Bolzano-Bozen, Faculty of Science and Technology, Bolzano (1) - Chiang Mai University, Department of Product Development Technology, Chiang Mai (2)

Keywords: PTR-QMS; GC-MS; coffee powder; aroma; volatile organic compounds

Introduction

Coffee is one of the most consumed beverages all over the world, as well as in Thailand especially in recent years. Two species of coffee are harvested in Thailand, *Coffea arabica* in the northern part and *Coffea canephora* var. *robusta* in the southern part. Arabica coffee plantation in the north of Thailand cover about 9000 hectares, distributed in eight provinces. Arabica coffee is mainly grown in the cooler highlands, wet processed in order to obtain high quality green bean and used in roasted coffee for domestic use[1].

Wet-processed beans are characterized by more fruity and acidic notes and less bitter and burnt attributes, compared with dry-processed beans. Moreover, the wet method allows the preservation of intrinsic properties of the beans, leading to coffees with a more desirable aroma[2].

Coffee aroma is the main responsible for organoleptic properties and therefore appeal of coffee. It is developed during the roasting process and depends on operative conditions such as temperature and time, but also on variety, origin and composition in terms of aroma precursors in green beans. In particular, caramelization of sugars leads to the production of furans (responsible for burnt and caramel notes), Maillard reaction of sugars develops pyrazines (characterized by roasty notes, mainly responsible for the coffee aroma), trigonellines are the precursors of pyridines (which confer smoky notes to the coffee) and lipids metabolism leads to the production of carbonyls [2, 3, 4, 5].

The effect of three different factors, such as altitude of harvesting (900 and 1500 meters above sea level), germination time (0, 24, 28 and 72 hours) and roasting level (light, medium and dark) on the volatile profile found in 24 different types of coffee were analyzed by Proton Transfer Reaction Quadropole-Mass Spectrometry (PTR-QMS) and by Gas Chromatography Mass Spectrometry (GC-MS). The best conditions to minimize undesirable odor and maximize the desirable ones were evaluated by Design Expert®.

Methods

Coffee beans (*C. arabica*) were harvested at 900 and 1500 meters above sea level in the Thai province of Chiang Rai in Doi Chang by Waweecoffee Co., Ltd. 1 Kg of green coffee beans, harvested at each altitude, was placed in boxes, soaked with 5 liters of water in a growth chamber at 28°C, in dark, and let germinate for 24, 48 and 72 hours. 300 gr of green coffee beans harvested at 900 and 1500 meters were not induced to germinate. Germinated and not germinated coffee beans were roasted at 200°C by Waweecoffee Co., Ltd. Three different roasting profile were achieved: light, medium and dark. Light-roasted coffee beans were roasted for 6 minutes, show 1-5% of weight loss and a light brown color; medium-roasted coffee beans were treated for 8 minutes, show 5-8% of weight loss and a medium brown color; dark-roasted coffee beans were roasted for 10 minutes, show 8-10% of weight loss and a dark brown color.

The headspace of coffee powders was measured by PTR-QMS. 35 mg of coffee powder were placed to equilibrate in 40-mL glass vials at 30°C for 15 min. Five replicates for each sample were analyzed. Four empty vials were analyzed and considered as blanks. The measurement order was randomized to avoid possible memory effect. The instrumental conditions were the following: drift voltage 600 V, drift temperature 75°C, drift pressure 2.10 mbar, affording an E/N value of 149 Td ($1 \text{ Td} = 10^{-17} \text{ cm}^2 \text{ V}^{-1} \text{ s}^{-1}$).

SPME analysis was carried out with a Gas Chromatography-Mass Spectrometer. 2 g of coffee powder were weighted in a 20-mL septum-sealed gas vial and sampled through the SPME fiber for 30 min at 60°C. Chromatographic conditions were the following: 40°C for 3 min, from 40°C to 200°C (0 min) at 3°C/min, then to 240°C (1.2 min) at 16°C/min.

Multivariate data analysis methods were applied to visualize data and the optimization of factors was performed using Design Expert®.

Results

2 g of coffee powder from eight different samples (L900, D900, L924, D924, L1500, D1500, L1524 and D1524) were analyzed by SPME-GC-MS. 27 compounds were tentatively identified and five of them were found only in the high-altitude samples. Acetic acid and 2-furanmethanol were high in L900 and L1500 samples and decreased with roasting and germination. 2-ethyl-6-methyl-pyrazine is found to be more abundant in low-altitude samples compared to the high-altitude ones and show no significant changes

related to roasting level and germination time. Dodecane undergoes a reduction in both low- and high-altitude samples after 24 hours of germination

The volatile profile of all 24 types of coffee was monitored by PTR-QMS. The 181 mass peaks obtained by each of these analyses were subjected to one-way ANOVA ($p < 0.01$); subsequently 89 masses were found at significantly different concentration between the three factors. Among these 89 masses, only the ones also found in GC-MS and tentatively identified (14 masses) were used in the optimization of the conditions. A 3-level factorial design was used to study the effect of altitude, germination time and roasting level in the composition of desirable and undesirable volatile compounds in coffee. The 14 masses found both in GC-MS and in PTR-MS were used in the model and the optimum conditions were found minimizing undesirable odor (m/z 61, 80, 99 and 101), and maximizing the desirable ones (m/z 95, 97, 109, 111, 117, 123, 131, 135, 137 and 148). The software proposed 23 solutions, with an index of desirability between 0.537 (best solution) and 0.256 (worst one). The two best solutions were the samples M924 (low-altitude, 24h-germination and a medium level of roasting) with an index of desirability of 0.537 and M1524 (high-altitude, 24h-germination and a medium level of roasting), with an index of desirability of 0.517.

Conclusions

Altitude of harvesting, germination time and roasting level have different impact in the volatile profile of coffee powder. Altitude of harvesting seems to be the factor contributing the less, while germination time and roasting level have a more important role in the development of aroma. The tentative identification of mass peaks identified by PTR-QMS and GC-MS allows the determination of desirable and undesirable odor and therefore the optimization of the three factors responsible for the volatile profile of coffee powder.

Novel Aspect

PTR-QMS and GC-MS has been used for the first to evaluate the effect of altitude of harvesting, germination time and roasting level in the development of aroma profile in coffee powder. Moreover, it is the first time in which data from PTR-QMS and GC-MS has been used coupled with Design Expert® to optimize growing conditions and treatments of coffee beans to obtain a more desirable aroma in coffee powder.

References

1. Angkasith P., AU Journal of Technology, 3, 117-123 (2002).
2. Lee L. W., Cheong M. W., Curran P., Yu B., Liu S. Q. Food Chemistry, 185, 182-191(2015).
3. Ruosi M. R., Cordero C., Cagliero C., Rubiolo P., Bicchi C., Sgorbini B., Liberto E. Journal of Agricultural and Food Chemistry, 60, 11283-11291(2012).
4. Yeretian C., Jorband A., Lindinger W. International Journal of Mass Spectrometry 223-224, 115-139 (2002).
5. Cheong W. M., Tong K. H., Ong J. J. M., Liu S. Q., Curran P., Yu B. Food Research International, 51, 388-396 (2013).

MP-124 / FEASIBILITY OF SELECTED ION FLOW TUBE MASS SPECTROMETRY (SIFT-MS) TO REAL-TIME MONITORING OF EXHALED BREATH VOLATILE ORGANIC COMPOUNDS (VOCS) IN CANCER PATIENTS.

Giuseppe Corona⁽¹⁾ - *Gianmaria Miolo*⁽²⁾ - *Andrea Caretta*⁽³⁾ - *Agostino Steffan*⁽¹⁾

National Cancer Institute Immunopathology and Cancer Biomarkers Aviano, Italy (1) - National Cancer Institute Medical and Preventive Oncology Aviano, Italy (2) - SRA Instruments Analytical Division Cernusco sul Naviglio, Italy (3)

Keywords: Selected Ion flow tube mass spectrometry, Breath, Real-time analysis, Volatile organic compounds, Cancer.

Introduction

The VOCs in the exhaled breath represent promising biomarkers able to give useful information on the physio-pathological status of the individual [1]. Real-time breath analysis of VOCs is an attractive non-invasive approach for cancer disease diagnosis and monitoring. Here we report the feasibility study of the use of selected ion flow tube mass spectrometry (SIFT-MS) to measure targeted breath VOCs suitable to monitoring cancer patients in a clinical setting.

Methods

The analysis of 32 targeted VOCs has been performed by Voice 200 Ultra™. The instrumental VOCs panel included: ammonia, acetone, isoprene, acetonitrile, dimethyl ether, aliphatic alcohols (n=5, C₁-C₅), cyanuric acid, low chain aliphatic acids and their corresponding aldehydes (n=12, C₁-C₆) including benzaldehyde, sulfur ethers (n=5) and phenol compounds (n=3). The study involved (n=12) healthy controls and (n=17) metastatic cancer patients.

Results

Twenty two (72%) out of 32 VOCs targeted VOCs panel were detectable by SIFT-MS under the analytical conditions used. As compared to healthy group cancer patients present significant lower breath level of ammonia (p<0.001), acetone (p<0.05), carbon disulfide (p<0.01) and phenol (p<0.05) evaluated by t-test. Within the series investigated was possible to recognize qualitative and quantitative individual VOCs signatures for both cancer patients and healthy individuals. Specific VOCs signature associated with the disease progression was observed in a case of a patient with sarcoma.

Conclusions

The use of SIFT-MS in clinical setting is a feasible and appropriate approach for the qualitative and quantitative analysis of a discrete number of exhaled VOCs. The study underlines the potential of the real-time VOCs profile analysis to identify specific patient pharmacological phenotype and supports the use of such breath metabolomics approach for longitudinal monitoring of cancer patients.

Novel Aspect

The SIFT-MS analytical instrumentation can be effectively applied for the real-time breath VOCs profile analysis in a clinical setting

References

1. A. Amann, W. Miekisch, J. Schubert, B. Buszewski, T. Ligor, T. Jezierski, J. Pleil, and T. Risby, *Annu. Rev. Anal. Chem.* 7. 455-482, 2014.

MP-125 / CLOSTRIDIUM TETANI (C. TETANI) CULTIVATION PROCESS MONITORING BY HEADSPACE SOLID-PHASE MICROEXTRACTION (HS-SPME) COMBINED WITH GC-MS

Ali Es-haghi⁽¹⁾ - *Masoud Ghader*⁽¹⁾ - *Nader Shokoufi*⁽²⁾ - *Kazem Kargosha*⁽²⁾

Razi Vaccine & Serum Research Institute, Agricultural Research, Education and Extension Organization (AREEO) Department of Physico Chemistry, Karaj (1) - Chemistry & Chemical Engineering Research Center of Iran, Faculty of Clean Technologies, Tehran (2)

Keywords: Process monitoring; volatile organic compounds; headspace solid-phase microextraction; GC-MS; C. tetani.

Introduction

Vaccines are one of the most important products of biological processes which are generally less controllable and more susceptible to variation [1]. HS-SPME-GC-MS method, allows determination of biogenic volatile organic compounds (BVOCs) emerged from microorganism [2]. The identified volatiles can

be considered as biomarkers for the identification and classification of organisms as well as for monitoring the bioreactor [3-5].

Methods

HS-SPME sampling was applied for extraction of VOCs, emerged from *C. tetani* in different culture media, followed by GC-MS analysis for identification and quantitation of BVOCs. Then, a portable dynamic air sampling (PDAS) device was constructed and used for sampling biomarker VOCs from the exhaust gas flow of bioreactor by the SPME fiber at the optimum condition. Other cultivation parameters such as pH and optical density were measured during the sampling time.

Results

About fifty volatile compounds released from *C. tetani* in different media were identified by the proposed method. The main detected compounds were sulfur compounds such as dimethyl trisulfide, dimethyl tetrasulfide, and butanethioic acid S-methyl ester which specified as a biomarker for monitoring bacterial growth in a bioreactor.

In the second step, sampling condition was optimized and the exhaust gas flow rate of 1mLs⁻¹ and sampling time of SPME fiber in PDAS of 20 minutes were obtained. Sampling of bioreactor was performed 6 times per day during the incubation period. The results showed good correlation between conventionally control parameters of bioreactor and quantitated VOCs.

Conclusions

Dynamic sampling from off-gas of bioreactor by SPME fibers was used for sampling bacterial VOCs from bioreactor and analysis was done by GC-MS. Reliability of proposed method was performed using pH and optical density of the *C. tetani* cultivation as the conventional monitoring parameters. results showed that, the proposed method can be used as a novel process monitoring method for microorganism cultivation.

Novel Aspect

HS-SPME followed by GC-MS method was developed for monitoring bacteria growth in bioreactor.

References

1. Streefland M., van de Waterbeemd B., Happé H., van der Pol L. A., Beuvery E. C., Tramper J., Martens D. E., Vaccine, 25, 2994–3000 (2007).
2. Zhang Z., Li G., Microchemica Journal, 95, 127-139 (2010).
3. Tait E., Perry J. D., Stanforth S. P., Dean J. R., Journal of Chromatographic Science, 52, 363-373(2014).
4. Liu Y., Miao Z., Guan W., Sun B., Molecules, 17, 3708–3722(2012).
5. C. de V. Petersen, H. C. Beck, F. R. Lauritsen, Biotechnology and Bioengineering, 85, 298–305, (2004).

MP-126 / SPME-GC/MOX SENSOR SYSTEM: A NEW METHOD FOR THE EVALUATION OF HUMAN SEMEN QUALITY

Valentina Longo⁽¹⁾ - **Angiola Forleo**⁽¹⁾ - **Sara Pinto Provenzano**⁽²⁾ - **Daniela Domenica Montagna**⁽²⁾ - **Lamberto Coppola**⁽²⁾ - **Vincenzo Zara**⁽³⁾ - **Alessandra Ferramosca**⁽³⁾ - **Pietro Siciliano**⁽¹⁾ - **Simonetta Capone**⁽¹⁾

National Research Council, Institute for Microelectronics and Microsystems, Lecce (1) - Biological Medical Center “Tecnomed”, Biological Medical Center “Tecnomed”, Nardò (2) - University of Salento, Department of Environmental and Biological Sciences and Technologies, Lecce (3)

Keywords: volatile compounds, SPME, GC-MOX sensor system, human semen

Introduction

The most used method for determining the volatile compounds of biological matrices is headspace solid-

phase microextraction gas chromatography (HS-SPME-GC). Usually, the detector is mass spectrometer (MS) [1]. In this work, we propose a new technique based on the application of a MOX sensor [2], as additional detector after GC separation, to evaluate metabolome profile of human semen applied to infertility study.

Methods

First, semen samples motility was assessed by standard semen analysis (SA). Next, volatile organic compounds (VOCs) of human semen were extracted by SPME technique and separated by gas chromatograph. The two detectors, MS and MOX sensor, work in parallel by a two-way splitter. The MOX sensor traces (resistance vs time) and physiological parameters of semen analysis were used for data analysis. Partial least-squares discriminant analysis (PLS-DA) was applied to the preprocessed sensorgrams to evaluate the ability of the GC/MOX sensor to discriminate sperm samples with different motility.

Results

Seminograms allowed us to classify the 57 samples in two motility groups: 16 asthenozoospermic samples with sperm progressive motility <32% and 41 normozoospermic ones. Forty-nine total VOCs were detected by GC/MOX sensor system.

Overlapping of a chromatogram and the corresponding sensorgram (as function of run time) for each sample group showed that background contamination, such as column bleed, is not perceived by VOC sensor.

2D- and 3D- PLS-DA results obtained from sensor data were compared with data obtained by seminograms, considering 2 groups (asthenozoospermic and high motility samples). 2D-PLS-DA based on SA showed an accuracy of 85% in the discrimination in two groups. In three dimensions, accuracy increases up to 92%. Using exclusively sensor resistance profiles, accuracy is of 77% (in 2D) and 78% (in 3D) for discrimination in asthenozoospermic and normozoospermic samples.

VOC sensor has different affinity to organic compound classes, having a higher reactivity with aldehydes, ketones and acetamides.

Conclusions

GC/MOX sensor system allowed us to develop a reproducible and novel method to detect human semen headspace VOCs. We have shown that the never-explored seminal VOC fingerprinting is a valid basis for discriminating samples with different motility. PLS-DA based on VOCs analysis by sensorgrams provided good classification between asthenozoospermic from normozoospermic samples, compared to PLS-DA based on sperm motility measured by seminogram.

Novel Aspect

Setup of a dual system of volatile compound detection that allow to outline resistance profile that can discriminate asthenozoospermic from normozoospermic samples.

References

1. García-Llatas G., Lagarda M. J., Romero F, AbelláP., Farré R, Food chemistry, 101, 1078-1086 (2006).
2. Longo V., Forleo A., Pinto Provenzano S., Montagna D. D., Coppola L., Zara V., Ferramosca A., Siciliano P., Capone S. Sensors. Proceeding of CNS2018. (2018), in press.

MP-127 / THE VOLATILE FAECAL METABOLOME IN CUSHING'S SYNDROME MOUSE MODEL OF TYPE 2 DIABETES

Célia Lourenço⁽¹⁾ - ***Michael Cauchi***⁽²⁾ - ***Liz Bentley***⁽³⁾ - ***Roger D. Cox***⁽³⁾ - ***Claire Turner***⁽⁴⁾

Cranfield University, School of Aerospace, Transport and Manufacturing, Cranfield, Bedfordshire (1) - University of Limerick, Department of Mathematics & Statistics, Limerick (2) - MRC Harwell, Harwell Science and Innovation Campus, Mammalian Genetics Unit, Oxfordshire (3) - The Open University, School of Life, Health & Chemical Sciences, Faculty of Science, Technology, Engineering and Mathematics, Milton Keynes (4)

Keywords: Volatile organic compounds, metabolomics, gut microbiome, SIFT-MS, TD-GC-MS.

Introduction

Increasing evidence suggests that the host-gut microbiota interactions and gut microbiota composition impacts strongly on the metabolic phenotypes of the host, although the underlying mechanisms are poorly understood [1]. Dysbiosis of the gut bacteria influences the production of the volatile fermentation products resulting on an altered metabolome and contributing to the onset of metabolic diseases [2].

Methods

This research aimed to investigate longitudinal changes in volatile organic compounds (VOCs) present in the faecal headspace of Cushing's syndrome mouse model of type 2 diabetes and once compared to their wild-type littermates. The volatile faecal metabolome was analysed using Selected Ion Flow Tube Mass Spectrometry (SIFT-MS) and Gas Chromatography-Mass Spectrometry (GC-MS) combined with thermal desorption (TD).

Results

Univariate and multivariate data analysis revealed statistically significant differences in volatile metabolites between wild-type littermates and diabetic animals. Cushing's mice developed obesity and showed significantly higher blood glucose concentrations and plasma insulin concentrations, supporting the hypothesis that Cushing's mice develop a diabetic phenotype related to a gut microbiota interaction.

Conclusions

The gut microbiota of Cushing's diabetic mice seem to have a different composition once compared to wild-type littermates, and the statistical analysis revealed statistically significant differences in volatile faecal metabolites. Overall, the predictive model provided information on whether the volatile metabolome is able to discriminate groups according to the phenotype.

Novel Aspect

The use of animal models can aid in the interpretation of metabolic pathways under controlled conditions. This work has never been done before.

References

1. Sommer F., Bäckhed F., Nat. Rev. Microbiol., 11, 227-238 (2013).
2. Carding S., Verbeke K., Vipond D. T. Corfe B. M., Owen L. J., Microb. Ecol. Health Dis., 26, 1-9 (2015).

MP-128 / INVESTIGATIONS OF VOCs PRESENT IN BREATH USING PTR-MS: PRELIMINARY STUDIES OF KETONES AND ANAESTHETICS

Michaela Malásková⁽¹⁾ - *David Olivenza-Leon*⁽²⁾ - *Felix Piel*⁽³⁾ - *Pawel Mochalski*⁽¹⁾ - *Martin Klieber*⁽¹⁾ - *Christopher Mayhew*⁽¹⁾

Institute for breath research, University of Innsbruck, Innsbruck (1) - School of Physics and Astronomy, University of Birmingham, Birmingham (2) - Ionicon Analytik Gesellschaft m. b. H., Eduard-Bodem-Gasse 3, Innsbruck (3)

Keywords: breath, volatiles, proton transfer reaction mass spectrometry

Introduction

The need for a non-invasive clinical diagnostics has grown over the past years. Breath research is gaining an important role in this field. More than 500 volatile organic compounds (VOCs) in breath can provide us with a wide variety of biomarkers. Most of these compounds are present in very low concentrations (sub parts per billion), therefore finding a sensitive instrument for analysing breath samples is crucial.

Methods

Proton transfer reaction mass spectrometry (PTR-MS) has proven its qualities in many research fields including breath research. (1,2) Its high sensitivity for VOCs, easy or usually no sample preparation and if used correctly, high selectivity, provide us with a strong analytical tool. Although PTR-MS usually uses H_3O^+ as a reagent ion (RI), it can also be used in modes to produce different RIs, e.g. O_2^+ and NO^+ , which react with VOCs of interest providing the reactions are exothermic. 3

Results

The investigation with the ketones is a major collaborative project involving (3) ESRs from 3 organisations. We measured and analysed fragmentation patterns while changing the reduced electric field for 21 compounds. (4) A fast GC was coupled to the PTR-MS in order to separate isobaric compounds, which can have often completely different origins in the human body and therefore can be found in the breath. Some of these substances were measured also with different RIs and different humidity conditions, to provide us with an even wider database to allow ease of identification under different conditions.

In addition to investigating endogenous breath compounds, we have begun to investigate important exogenous compounds in breath. To date we have focused on the anaesthetics, isoflurane, sevoflurane, desflurane and enflurane. Initial studies have been undertaken to determine the product ions following reaction with H_3O^+ , in order to determine which ions should be monitored. Detailed studies of product ions as a function of reduced electric field will be presented.

Conclusions

We demonstrate that changing the reduced electric field aids in the identification of VOCs. This pre-study has generated a useful database of fragment product ions resulting from the ketones and anaesthetics for various reagent ions and under various conditions, which will be presented.

Novel Aspect

In our work we are aiming to create a database of VOCs present in breath and their reactions with RIs, the product ions of which could serve as biomarkers for diseases.

References

1. Ilis A. M., Mayhew C. A., John Wiley & Sons Ltd., UK, 2014.
2. Blake R. S., Monks P. S., Ellis A. M., Chem. Rev., 109, 861–896, 2009.
3. Spanel P., Ji Y., Smith D., International Journal of Mass Spectrometry and Ion Processes, 165/166, 25-37, 1997.
4. Spanel P., Smith D., International Journal of Mass Spectrometry, 280, 128–135, 2009.

MP-129 / AN INNOVATIVE ANALYSIS METHOD FOR FRAGRANCE OF SAKE USING AN INTRODUCING DEVICE FOR VOLATILE COMPOUNDS COMBINED WITH DART-MS

Chikako Takei⁽¹⁾ - *Yukiko Otake*⁽¹⁾ - *Kenichi Yoshizawa*⁽¹⁾ - *Haruka Nishimoto*⁽²⁾ - *Takahiro Akashi*⁽²⁾

BioChromato, Application support, Fujisawa (1) - Hakutsuru Sake Brewing Co., Ltd., Research & Development, Kobe (2)

Keywords: Sake, Fragrance, DART-MS

Introduction

Sake is a Japanese national liquor, and many kinds are brewed for different flavor notes such as fruits, flowers, herbs, and spices. Additionally, sake can be served in a variety of vessels (cups): a choko, a Masu, a sakazuki, and a wineglass, allowing one to enjoy different scents, depending on the cup. The purpose of this work is to visualize the differences of sake fragrance depending on cups.

Methods

Hakutsuru Daiginjo (sake) and cups (Choko, Daichoko, wineglass, and Champagne glass) were used. Analytical system was composed of compact (Bruker) mass spec equipped with DART ion source (IonSense) and Volatimeship introducing device for volatile compounds (BioChromato). Data acquisition via the mass spectrometer were began, the cup added in sake were placed at the introducing tube of the Volatimeship, where the position of human nose, when the human is drinking.

Results

Ethyl caproate (*i.e.*, ethyl hexanoate) provides one of the characteristic fragrances of sake. As the cup of sake was placed at the introducing tube of the Volatimeship, ethyl caproate was immediately detected by mass spectrometer. The volatilization behavior of ethyl caproate for each cup were measured. The wineglass not only showed higher intensity at the moment of placing the cup, but it also showed better sustainability. These analysis results supported the results of organoleptic evaluation. Moreover, the detected compounds from sake was different depend on cups. As for wineglass and champagne glass, there were some compounds that were not detected from Choko. These analysis results also supported the results of organoleptic evaluation, too.

Conclusions

An introducing device for volatile compounds combined with DART-MS enabled continuous detection of the change in fragrance intensity directly from the sake vessel. Therefore, this method can be useful for objective evaluation of the duration and change of scent. Additionally, this method can be useful for visualizing the sensory evaluation, because it enabled visualization of the varying intensity, or quality, of the fragrance of sake depend on cups.

Novel Aspect

An introducing device for volatile compounds combined with DART-MS enables objectively evaluation the sensual sensation of fragrance.

MP-130 / EXPLORING INTER-INDIVIDUAL DIFFERENCES IN IN-VIVO FLAVOUR RELEASE BY SIFT-MS

Michele Pedrotti^(1,2) - ***Leonardo Menghi***⁽¹⁾ - ***Iuliia Khomenko***^(1,3) - ***Andrea Caretta***⁽⁴⁾ - ***Eugenio Aprea***⁽¹⁾ - ***Annachiara Cavazzana***⁽⁵⁾ - ***Thomas Hummel***⁽⁵⁾ - ***Flavia Gasperi***⁽¹⁾ - ***Franco Biasioli***⁽¹⁾

Fondazione Edmund Mach, Research and Innovation Center, Department of Food Quality and Nutrition – via E. Mach 1, 38010 San Michele all'Adige (TN, Italy) (1) - Food Quality & Design Group, Wageningen University & Research, WG NL-6708, Wageningen, (Netherlands) (2) - Institute for Ion Physics and Applied Physics, University of Innsbruck, Technikerstr. 25, Innsbruck (Austria) (3) - SRA Instruments S. p. A., via alla Castellana 3 - 20063 Cernusco sul Naviglio (MI, Italy) (4) - Interdisciplinary Center Smell & Taste Department of Otorhinolaryngology TU Dresden, Medical Faculty Carl Gustav Carus, Dresden (Germany) (5)

Keywords: nose-space analysis, in vivo flavor analysis, SIFT-MS, population screening

Introduction

Selected-Ion Flow-Tube Mass Spectrometry (SIFT-MS) allows on line accurate quantification of several volatile organic compounds down to the ppb range [1]. The technique has been applied to different cases of breath analysis [2]–[5] but to the authors knowledge it has not been applied to monitor in nose flavor release during food consumption for a large population.

In this work we describe a SIFT-MS based approach to measure in-vivo nose-space VOCs concentration during consumption of strawberry flavored candies on a large healthy population with the aim to investigate the inter-individual variability in flavor release: nose space analysis might disclose new information regarding the interaction between products and consumers [6].

Methods

In this study the effect of gender, age and BMI on flavor release was investigated. Ninety-two volunteers (60.9 % female) between 22 to 68 years old (mean= 40.9; SD= 11.9) participated in the study. Each subject, after filling in a socio-demographic questionnaire, consumed a strawberry candy following a specific bite-based procedure supported by a video tool. Simultaneously nose-space analysis with a SIFT-MS machine was carried out (SYFT VOICE 200 ultra, Syft Ltd, New Zealand). Each subject performed at least 3 replicates of the same sample.

Results

With the set up used it was possible to monitor in real time 7 different aroma compounds released by the candy (ethyl maltol, 3-hexen-1-ol, ethyl 2-methylbutanoate, (Z)-3-hexenyl acetate, ethyl butanoate, ethyl hexanoate, 2-methylbutanoic acid) and acetone as marker of the participant breathing cycle. Principal Component Analysis (PCA) showed a possible effect of age and BMI on the flavor release signals. Both age and BMI were found to have a slightly negative correlation with in-nose aroma release for three of the monitored compounds (Pearson correlation, $p < 0.05$). Significant differences were found in hexenyl acetate and acetone levels (Welch t-test $p < 0.01$) between males and females: in both cases females presented higher in-nose levels. The study is part of a broader investigation that aims to evaluate the Italian olfactory function (COLFIT PROJECT) and here some preliminary results of the explorative data analysis are presented.

Conclusions

Real time mass spectrometry analysis by SIFT-MS can detect and monitor in-vivo volatile organic compounds release during candies consumption. For the first time the analysis was applied to a large sample of the population with the aim to characterize how physiological factors like age, gender, BMI and eating behaviors affect in nose aroma release. The results suggest an effect of age on aroma release and an effect of BMI that might be used to study some aspects underlying the retronasal perception mechanisms in more detail.

Novel Aspect

Real time direct injection mass spectrometric analysis through SIFT-MS is a powerful tool to realize population screening and investigate flavor release mechanisms.

References

1. D. Smith, P. Španěl, "Selected ion flow tube mass spectrometry (SIFT-MS) for on-line trace gas analysis," *Mass Spectrom. Rev.*, vol. 24, no. 5, pp. 661–700, 2005.
2. P. Španěl, D. Smith, "Progress in SIFT-MS: Breath analysis and other applications," *Mass Spectrom. Rev.*, vol. 30, no. 2, pp. 236–267, 2011.
3. B. Enderby, W. Lenney, M. Brady, C. Emmett, P. Španěl, D. Smith, "Concentrations of some metabolites in the breath of healthy children aged 7-18 years measured using selected ion flow tube mass spectrometry (SIFT-MS)," *J. Breath Res.*, vol. 3, no. 3, 2009.
4. P. Španěl, D. Smith, "Quantification of trace levels of the potential cancer biomarkers formaldehyde, acetaldehyde and propanol in breath by SIFT-MS," *J. Breath Res.*, vol. 2, no. 4, 2008.
5. R. N. Bloor, P. Španěl, D. Smith, "Quantification of breath carbon disulphide and acetone following a single dose of disulfiram (Antabuse) using selected ion flow tube mass spectrometry (SIFT-MS)," *Addict. Biol.*, vol. 11, no. 2, pp. 163–169, 2006.
6. A. Romano et al., "Nosespace analysis by PTR-ToF-MS for the characterization of food and tasters: The case study of coffee," *Int. J. Mass Spectrom.*, vol. 365–366, pp. 20–27, 2014.

MPS-S06 - Clinical chemistry

Chairs: Hugo Rocha, Giancarlo la Marca

MP-131 / IMPACT OF CYP2D6 GENETIC POLYMORPHISM ON THE PHARMACOKINETICS OF TRAMADOL AND ITS THREE MAIN METABOLITES IN KOREANS

Jung-Woo Bae - Sooyeon Lee - Kyung-Soo Chun - Min Je Choi - Seong-Kuk Hong

Keimyung University, College of Pharmacy, Daegu

Keywords: Tramadol, CYP2D6, LC-MS/MS

Introduction

Tramadol is a centrally acting analgesic used in the treatment of moderate to severe pain. Tramadol is mainly metabolized by the CYP2D6, CYP2B6 and CYP3A4 to *O*-desmethyltramadol (ODMT), *N*-desmethyltramadol (NDMT) and *N,O*-didesmethyltramadol (NODMT) [1]. We investigated the effects of CYP2D6 genetic polymorphism on the pharmacokinetics of tramadol and its three main metabolites.

Methods

A multiple oral dose of tramadol was administered to twenty-two subjects with CYP2D6EM (n=14) and CYP2D6IM (n=8). Blood samples were collected up to 72 hr after drug intake, and plasma concentrations of tramadol and its metabolites (ODMT, NDMT and NODMT) were determined by validated HPLC-MS/MS method.

Results

AUC of tramadol ($p<0.001$), NDMT ($p<0.05$) and NODMT ($p<0.001$) in CYP2D6IM group was significantly higher than those in CYP2D6EM group. The AUC ratios of ODMT to tramadol in the CYP2D6IM group was decreased to 44.9% of that in the CYP2D6EM group ($p<0.001$). The AUC ratios of NDMT to tramadol in the CYP2D6IM group was increased to 148.8% of that in the CYP2D6EM group (N. S.).

Conclusions

The CYP2D6 genotype was found to play a significant role in the metabolic conversion of tramadol into ODMT and, as a compensatory consequence, impacted the conversion of tramadol into NDMT.

Novel Aspect

CYP2D6*10 allele is found to affect the pharmacokinetics of tramadol and its phase I metabolites in Koreans.

References

1. Gong L., Stamer UM., Tzvetkov MV., Altman RB., Klein TE., Pharmacogenetics and Genomics, 24, 374-380 (2014).

MP-132 / INDIVIDUAL AND CYCLIC ESTROGENIC PROFILE IN WOMAN: STRUCTURE AND VARIABILITY OF THE DATA

Cristina Bozzolino - Sara Vaglio - Eleonora Amante - Eugenio Alladio - Enrico Gerace - Marco Vincenti

Centro Regionale Antidoping e di Tossicologia "A. Bertinaria", University of Turin, Department of Chemistry, Turin

Keywords: GC-MS, estrogens, natural rhythmic fluctuations, multivariate data analysis, female reproductive cancers

Introduction

In premenopausal women the role of estrogens is not thoroughly understood due to the large variations in estrogenic concentration throughout the menstrual cycle [1-6]. Lifestyle factors, such as smoking, physical activity and alcohol consumption may affect menstrual cycle outcomes [7]. Vice versa, the variations could influence fertility [8], the sexual behavior [9] and some diseases, such as schizophrenia [10], stress [11], insomnia [12] and anxiety [13].

Methods

A Gas Chromatography coupled to Mass Spectrometry (GC-MS) method capable of detecting 15 estrogens (17 α -estradiol, 17 β -estradiol, estrone, estriol and their principal metabolites) in urine was developed following the experimental design approach and fully-validated in accordance with the 17025:2005 principles. Multivariate data analysis was performed applying PARAFAC, Principal Component Analysis (PCA) and Linear Discriminant Analysis (LDA) to the data.

Results

The following validation parameters were successfully investigated: linearity range, selectivity, specificity, LOD, LOQ, intra-assay precision and accuracy, repeatability, matrix effect, extraction recovery and carry-over.

Urine samples derived from 9 childbearing female volunteers (age between 25-35) who collected their urine every day for 28 days, that is a complete menstrual cycle. The data were evaluated by PARAFAC and the natural fluctuation of the concentration of each estrogen throughout the 28-day menstrual cycle was extract. Then, PCA was carried out on the data and the Scores Plot showed the occurrence of three different clusters corresponding to the three phases of the menstrual cycle: follicular phase, ovulation and luteal phase. Plotting PC1 and PC2 as a function of the menstrual cycle day in two separate graphs, the transitions between the phases turned visible by identifying the starting point of both the ovulation and the luteal phase. Lastly, a preliminary LDA model was built.

Conclusions

The application of a GC-MS method coupled to multivariate statistical techniques demonstrated great potentialities in the study of the natural rhythmic fluctuations of estrogens throughout the menstrual cycle. Therefore, the next step will be the application of this innovative approach to examine in depth the relationship between the endogenous estrogens and the increased risk of female reproductive cancers in postmenopausal women.

Novel Aspect

For the first time, not only the variation of the parent estrogens was studied, but also of their principal metabolites.

References

1. Brown S. B., Hankinson S. E., *Steroids*, 99, 8-10 (2015).
2. Ziegler R. G., Fuhrman B. J., Moore S. C., Matthews C. E., *Steroids*, 99, 67-75 (2015).
3. Mueck A. O., Seeger H., Lippert T. H., *Maturitas*, 43, 1-10 (2002).
4. De Almeida Chuffa L. G., Lupi-Júnior L. A., Balandis Costa A., de Arruda Amorim J. P., Ferreira Seiva F. R., *Steroids*, 118, 93-108 (2017).
5. Samavat H., Kurzer M. S., *Cancer Letters*, 356, 231-243 (2015).
6. Folkerd E., Dowsett M., *The Breast*, 22, S38-S43 (2013).
7. Liu Y., Gold E. B., Lasley B. L., Johnson W. O., *American Journal of Epidemiology*, 160, 131-140 (2004).
8. Barrett E. S., Parlett L. E., Windham G. C., Swan S. H., *Fertility and Sterility*, 101, 1773-1780 (2014).
9. Regan P. C., *The Canadian Journal of Human Sexuality*, 5, 145-156 (1996).

10. Bergemann N., Parzer P., Runnebaum B., Resch F., Mundt C., *Psychological Medicine*, 37, 1427–1436 (2007).
11. Roney J. R., Simmons Z. L., *Adaptive Human Behavior and Physiology*, 1, 30–40 (2015).
12. Driver H. S., Werth E., Dijk D. J., Borbély A. A., *Sleep Medicine Clinics*, 3, 1–11 (2008).
13. Ellison P. T., Lipson S. F., Jasienska G., Ellison P. L., *American Journal of Physical Anthropology*, 134, 513–519 (2007).

MP-133 / PHARMACOKINETIC OF MYRIOCIN IN RABBIT'S EYES

Michele Dei Cas⁽¹⁾ - ***Federico Maria Rubino***⁽¹⁾ - ***Jessica Rizzo***⁽¹⁾ - ***Signorelli Paola***⁽¹⁾ - ***Enrica Stretto***⁽²⁾ - ***Chiara Platania***⁽³⁾ - ***Claudio Bucolo***⁽³⁾ - ***Rosario Pignatello***⁽⁴⁾ - ***Riccardo Ghidoni***⁽¹⁾ - ***Rita Clara Paroni***⁽¹⁾

University of Milan, Department of Health Sciences, Milan (1) - CNR Neuroscience Institute, CNR Neuroscience Institute, Pisa (2) - University of Catania, Department of Biomedical and Biotechnological Sciences, Catania (3) - University of Catania, Department of Pharmaceutical Sciences, Catania (4)

Keywords: myriocin, pharmacokinetic, drug delivery, eye, qTrap

Introduction

Myriocin (Myr) is a suicide inactivator of ceramide synthesis with a complex lipid multifunctional structure. Its biological activity is exerted at very low doses, and thus highly performing quantitative method are needed [1]. The pharmacological development of Myr to modulate ceramide levels also requires currently unavailable ADME information in healthy and pathological animal models.

Methods

A 3200 qTRAP system was operated under manufacturer's instructions to perform extensive tandem MS studies. LC and MS were optimized for the separation, detection, identification and quantification of Myrin aqueous humour, vitreous humor and retina of rabbit after different treatments. MRM scan mode was used as protocol, to tackle the ambiguities of some samples, EPI scan, MS3 and the use of HR-MS on a 6600 qTOF was also crucial.

Results

Tissue levels of Myrare related to the dose, pharmaceutical formulation, administration route, and administration schedule. Four type of eye drop formulations were evaluated in rabbits: Myr in buffer solution, Myr in two different types of liposomes (neutral and cationic), Myr encapsulated in a novel type of solid nanolipid particles. In particular, solid nanolipid particles was found to be extremely efficacious for the delivery of Myrin the different areas of the eye. The slow release formulation produced a Tmax between 180-240 minutes, with a maximum concentration respectively of 10 µg/mL in vitreous and 1.7 µg/g in retina. Biological activity of Myr was always confirmed in rabbit's retina by a significant decrease in ceramides levels, compared to the physiological level established from untreated animals. Some untreated animals demonstrated the presence of a peak with spectroscopic and chromatographic characteristics of Myr. EPI scan, MS3, and the use of a HR-MS as protocol confirmed an upstream contamination.

Conclusions

A new pharmaceutical delivery system, such as solid lipid nanoparticles, allows Myr to accumulate into the retina by a non-invasive eye-drops administration, and to modulate sphingolipid levels. These results would be promising to study its efficacy in mouse pathological model of retinitis pigmentosa. MRM measurement needs alternative and/or complementary techniques to confirm eventual false-positives samples.

Novel Aspect

The use of HR-MS as protocol should be quite advantageous for predicting new metabolites and for

removing unwanted matrix interferences.

References

1. Campisi G.M., Signorelli P., Rizzo J., Ghilardi C., Antognetti J., Caretti A., Lazarević J.S., Strettoi E., Novelli E., Ghidoni R., Rubino F.M., Paroni R. *Biomed Chromatography*; 31:e4026. <https://doi.org/10.1002/bmc.4026> (2017)

MP-134 / COMPARISON OF ESI, APCI AND APPI IONIZATION TECHNIQUES FOR TESTOSTERONE (T), DIHYDROTESTOSTERONE (DHT) AND ESTRADIOL (E2) LCMS/MS ASSAY

Reena Desai⁽¹⁾ - **Pekka Keski-Rahkonen**^(1,3) - **D. Tim Harwood**⁽²⁾ - **David J. Handelsman**⁽¹⁾

ANZAC Research Institute, University of Sydney, Sydney, NSW 2139, Australia (1) - Cawthron Institute, Nelson 7010, New Zealand (2) - Current address: International Agency for Research on Cancer, F-69372 Lyon, France (3)

Keywords: Comparing ionization techniques for steroids

Introduction

High sensitivity measurement at low concentrations of non-derivatized testosterone (T), dihydrotestosterone (DHT) and estradiol (E2) requires optimizing ionization technique but few matrix-specific comparisons are available. We compared electrospray ionization (ESI), atmospheric pressure chemical ionization (APCI) and atmospheric pressure photoionization (APPI) for sensitivity and matrix interference in neat solutions and human serum.

Methods

LC instrumentation and conditions

A Shimadzu Prominence HPLC system was used and consisted of a DGU-20A5 degasser, three LC-20 AD pumps, a SIL-20AT HC autosampler and a CTO-20A column oven (Shimadzu Scientific Instruments, Kyoto, Japan). The column was an Agilent Poroshell 120 SB-C18 (2.1 × 50 mm; 2.7 μm) with 0.3 μm in-line filter (Agilent Technologies). For the determination of LLOQs with standard solutions, column temperature was 40°C, injection volume 20 μl, and mobile phase of 70% methanol/water, containing 0.1% (v/v) of formic acid for positive ESI and 2.5 mM ammonium hydroxide for negative ESI. APCI and APPI were used without mobile phase additions. For serum and tissue samples, a following gradient was employed (A: water, B: methanol): 0–1 min: 20% B, 1–8 min: 20 > 70% B, 8–10 min: 70% B, 10–12 min: 100% B, 12–16 min: 20% B. Flow rate was 0.3 ml/min, injection volume 40 μl and same mobile phase additions were used as for the isocratic separations. A flow-line selection valve was used to divert the eluent to waste during 0–7 min and 11–16 min. For all experiments with APPI, toluene was used as a dopant solvent, delivered with an LC-20 AD pump with a flow rate 10% of that of the mobile phase.

MS instrumentation and conditions

An API 5000 triple quadrupole MS was used with Turbo V (ESI, APCI) and PhotoSpray (APPI) ion sources (AB Sciex, Toronto, Canada). The APPI lamp was a DC-driven 10.0 eV krypton discharge lamp model PKS100 from Heraeus (Hanau, Germany). Nitrogen was employed as ion source, curtain and collision gas (Peak NM20ZL, Peak Scientific, Inchinnan, UK). For the precursor ion selection and optimization of the MS parameters, 5 μM solution of E2 in methanol was infused post-column into the mobile phase (70% methanol/water).

The E2 infusion flow rate was kept at 5% of the mobile phase flow. Precursor and product ion selection was based on ion intensity. After establishing the MRM conditions, further optimization of the ion source parameters and the mobile phase flow rate was made with on-column injections, which enabled signal-to-noise (S/N) calculations. MS parameters not related to ionization were optimized and kept identical for each MRM transition. Both quadrupoles were set at unit resolution, cycle time was 250 ms and no signal

thresholding was applied. Data acquisition and processing were performed with AB Sciex Analyst 1.5.2 software.

We used a Shimadzu Nexera UHPLC and SCIEX API 5000 triple quad with Phenomenex Kinetex XB C18 (50mm x 2.1mm x 1.7 μ m) column at 40 C, injection volume 50 μ l, flow rate was 0.5 ml/min, mobile phase of methanol-water. For ESI we used modifiers of 0.1% formic acid for positive and 2.5 mM ammonium hydroxide for negative polarity. APCI and APPI were used without mobile phase additions but toluene as a dopant for APPI.

Results

Direct infusions of 5 μ g/ml of T, DHT and E2 carried out to optimize the best ionization and MRM transitions in both positive and negative polarity for E2, whereas only positive mode for T and DHT with ESI, APCI and APPI. With ESI, peak responses increased with 0.1 % formic acid in positive mode and 2.5mM ammonium hydroxide in negative mode for E2 as well for DHT in positive ESI. APCI and APPI both resulted in water loss for main precursor ion for DHT and E2 with the highest intensity, although no modifiers were required for high signal intensity, toluene improved ionization with APPI. T had similar responses for all measurement. Regardless of ionization, methanol in the mobile phase enhanced ionization more than acetonitrile. After establishing the optimized MS/MS conditions serum extracts and neat samples were analyzed, to define LLOQs with a precision <20%, and accuracy of 80–120% and matrix effect. APPI had strongest S/N ratio and better LLOQs for all analytes compared to APCI and ESI.

Conclusions

Ionization technique and ion polarity both contribute to sensitivity and susceptibility to matrix interference. With serum extracts, ESI suffered from severe ion suppression compared to APPI and APCI. We conclude that irrespective of ionization, negative polarity was best for E2 whereas, for T and DHT positive polarity is best. APPI has the highest potential to produce high S/N ratio with least ion suppression.

To maximize the LC–MS assay sensitivity for E2 analysis, various approaches have been described. Sample preparation and LC separation of steroids in biofluids such as serum or plasma are reasonably well established, but the most efficient ionization technique for E2 is still to be determined. E2 can be ionized with all the common commercially available ionization APPI is much less prone to ion suppression than either APCI or ESI. Evidence indicates that ESI is especially susceptible to ion suppression effects. The high linear dynamic range and the applicability to a broad range of polarity makes APCI a good choice for steroids, especially when the APPI source is not available.

Novel Aspect

This study comparing ESI, APCI and APPI for the analysis of T, DHT and E2 in serum provides insight on achieving high sensitivity and choice for optimal ionization to avoid matrix interference.

Accurate measurement of estradiol (E2) is important in clinical diagnostics and research. High sensitivity methods are critical for specimens with E2 concentrations at low picomolar levels, such as serum of men, postmenopausal women and children. Achieving the required assay performance with LC–MS is challenging due to the non-polar structure and low proton affinity of E2.

Previous studies suggest that ionization has a major role for the performance of E2 measurement, but comparisons of different ionization techniques for the analysis of clinical samples are not available. In this study, female serum and endometrium tissue samples were used to compare electrospray ionization (ESI), atmospheric pressure chemical ionization (APCI) and atmospheric pressure photoionization (APPI) in both polarities. APPI was found to have the most potential for E2 analysis, with a quantification limit of 1 fmol on-column. APCI and ESI could be employed in negative polarity, although being slightly less sensitive than APPI. In the presence of biological background, ESI was found to be highly susceptible to ion suppression, while APCI and APPI were largely unaffected by the sample matrix. Irrespective of the ionization technique, background interferences were observed when using the multiple reaction monitoring transitions commonly employed for E2 (m/z 271 > 159; m/z 255 > 145). These unidentified interferences were most severe in serum samples, varied in intensity between ionization techniques and required efficient chromatographic separation in order to achieve specificity for E2.

References

1. Keski-Rahkonen P, Huhtinen K, Desai R, Harwood DT, Handelsman DJ, Poutanen M, Auriola S , LC-MS analysis of estradiol in human serum and endometrial tissue: Comparison of electrospray ionization, atmospheric pressure chemical ionization and atmospheric pressure photoionization. *J Mass Spectrom*, 2013 ,48:1050-1058
2. Harwood DT, Handelsman DJ, Development and validation of a sensitive liquid chromatography-tandem mass spectrometry assay to simultaneously measure androgens and estrogens in serum without derivatization, 2009, *Clin Chim Acta* 409:78-84

MP-135 / HIGHLY SELECTIVE AND FAST IDENTIFICATION OF BRUCELLA BY IMMUNOAFFINITY- MALDI-TOF MASS SPECTROMETRY

Hamideh Hamidi⁽¹⁾ - ***Mahshid Zarrineh***⁽¹⁾ - ***Ilnaz Soleimani***⁽¹⁾ - ***Ramin Bagheri Nejad***⁽²⁾ - ***Ali Es-haghi***⁽³⁾ - ***Alireza Ghassempour***⁽¹⁾

Medicinal Plants and Drugs Research Institute, Shahid Beheshti University, Tehran (1) - Razi Vaccine & Serum Research Institute (RVSRI), Agricultural Research, Education and Extension Organization (AREEO), Brucellosis Department, Karaj, Alborz (2) - Razi Vaccine & Serum Research Institute (RVSRI), Agricultural Research, Education and Extension Organization (AREEO), Department of Physico Chemistry, Karaj, Alborz (3)

Keywords: MALDI-TOF; Mass spectrometry; Brucella; Immunoaffinity

Introduction

The genus *Brucella* contains highly infectious species classified as biological threat agents [1]. Brucellosis is endemic in parts of Iran, and the annual incidence of the human and animal brucellosis has been reported [2]. The timely identification of the microorganism is necessary not only to treat patients effectively but also to establish outbreak management, source tracing, and threat analyses. Matrix-assisted laser desorption ionization–time of flight mass spectrometry (MALDI-TOFMS) is a rapid method used to analyze biological differences in microorganisms [3]. In this work, we have developed a highly selective method to quickly and accurately identify bacteria in milk/ blood samples combining MALDI-TOFMS and immunoaffinity enrichment/separation.

Methods

Synthesized monodispersed magnetic particles are activated and modified with protein A for oriented immobilization of monoclonal antibodies. 5 mL of human whole blood is spiked with *B. melitensis* cells. The spiked blood is injected into a blood culture bottle and cultured at 35°C in the presence of 10% CO₂. Antibodies-modified magnetic beads are used for extraction of *B. melitensis* from artificially inoculated milk and human blood samples. After inactivation of enriched bacteria, the ethanol/formic acid extraction procedure is applied and analyzed by MALDI-TOF MS.

Results

MALDI-TOF MS results for the identification of *Brucella* strains have been investigated after enrichment by the immunoaffinity technique based on immobilization of antibody on the magnetic beads. Initial spectra from *B. melitensis* biovar 16M have been obtained for pattern matching-based bacteria identification and subsequently the experimental spectrum from a blood/milk sample has been compared with the reference one. Because of the method's high sensitivity, the blood culture time required for diagnosis can be greatly reduced. As a proof, whole blood/ milk spiked with a low initial concentration (10² or 10³ cells per mL) of bacteria is cultured, then the antibody coated beads are added to these prepared samples after different blood culturing times. Finally, the bacteria attached on the surface of beads have been eluted and analyzed by the MALDI-TOFMS. Therefore, entire diagnosis process (treatment and mass-spec analysis of *Brucellae* from blood/milk media) can be accurately accomplished less than 48 hours by using this newly developed

method.

Conclusions

We have developed an immunoaffinity-MALDI-TOF MS method possessing a simple pretreatment and high sensitivity for rapid detection of *Brucella* spp. in humans and animal clinical specimens. This method indeed facilitates early accurate diagnosis of the disease which is required for further well-timed actions including effective anti-bacterial therapy or implementing control measures.

Novel Aspect

Immunoaffinity-MALDI-TOFMS method has been developed for the sensitive and rapid identification of *Brucella* strains.

References

1. Lista F. et al., BMC Microbiology, 11, 2011, 1471-2180 (2011).
2. Golshani M., Buozari S., Iran Biomed J., 21, 349-359 (2017).
3. Angeletti S., Journal of Microbiological Methods, 138, 20-29 (2017).

MP-136 / ANALYSIS OF PROTEIN BIOMARKERS FROM LOW SAMPLE VOLUMES

Trine Grønhaug Halvorsen - Cecilie Rosting - Maren C. S. Levernæs - Leon Reubsæet

School of Pharmacy, University of Oslo, Oslo

Keywords: protein analysis - antibodies – LC-MS/MS – affinity clean-up – biomarker analysis

Introduction

Targeted protein determination using liquid chromatography tandem mass spectrometry (LC-MS/MS) offers a more reliable alternative to traditional immunometric assays. The interest in applying LC-MS/MS in protein biomarker analysis is hence increasing. For low volume samples (<50 µL) a highly specific sample preparation step normally is necessary to detect low abundance biomarkers (low pM-level).

Methods

Two different approaches are discussed: 1) NanoLC-MS/MS analysis of human chorionic gonadotropin (hCG) from dried matrix samples (blood, serum, plasma and urine) by affinity sample clean-up of whole protein followed by tryptic digestion. 2) NanoLC-MS/MS analysis of progastrin releasing peptide (ProGRP) from 50 µL of serum using affinity extraction of proteolytic epitope peptide after digestion.

Results

Ovarian and testicular cancer marker hCG was detected from four different biological matrices (whole blood, serum, plasma and urine) from 15 µL of dried matrix spots at its upper reference level (low pM-level). The detection limits (S/N=3) ranged from 14.5 pM (whole blood) to 30.5 pM (urine). The method was briefly evaluated for dried serum samples providing satisfactory results (RSD(%)=13-29% and accuracy (97-106%)) [1].

The small-cell lung cancer marker ProGRP was also determined at its upper reference level (low pM-level). Detection limits (S/N=3) of 3.9 pM from 50 µL of serum was observed when using monoclonal antibodies developed against intact proteins to capture a proteotypic epitope peptide. Compared to intact protein extraction, this concept, which we call epitope fishing, provided cleaner extracts.

Conclusions

NanoLC and state-of-the-art triple quadrupole MS, and affinity clean-up is a powerful combination for determination of low abundance protein biomarkers available in low volumes. Both affinity clean-up of

whole protein (hCG in 15 μ L of dried matrix spot) and of proteotypic epitope peptide (ProGRP in 50 μ L serum) have been successfully enabled at low pM-level.

Novel Aspect

Determination of low abundance protein biomarkers is enabled from small volumes (< 50 μ L) by affinity clean-up and state-of-the art nanoLC-MS/MS.

References

1. Rosting C, Tran EV, Gjelstad A, Halvorsen TG. Determination of the low-abundant protein biomarker hCG from dried matrix spots using immunocapture and nano liquid chromatography mass spectrometry, J. Chromatogr. B 1077–1078 (2018) 44-51.

MP-137 / EARLY AND NON-INVASIVE MASS SPECTROMETRY TOOL FOR DIAGNOSING ASPERGILLOSIS

Radim Dobias⁽¹⁾ - **Anton Skriba**⁽²⁾ - **Tomas Pluhacek**⁽²⁾ - **Milos Petrik**⁽³⁾ - **Pavlina Lyskova**⁽⁴⁾ - **Vladimír Havlíček**⁽²⁾

Institute of Public Health, Partyzanske nam. 2633/7, Ostrava (1) - Institute of Microbiology, Academy of Sciences of the Czech Republic, Videnska 1083, Prague 4 (2) - Institute of Molecular and Translational Medicine, Palacky University, Olomouc (3) - Public Health Institute, Department of Parasitology, Mycology and Mycobacteriology, Usti nad Labem (4)

Keywords: aspergillosis; non-invasive; diagnosis; siderophore; imaging

Introduction

Aspergillus fumigatus is a ubiquitous saprophytic airborne fungus responsible for more than one million deaths every year. Its siderophores represent important virulence factors contributing to microbiome-metabolome dialog in a host. From diagnostic point of view, the application of *Aspergillus* secondary metabolites is promising due to non-invasiveness and high speed, sensitivity and potential for standardization.

Methods

Using a model of experimental aspergillosis in immunocompromised Lewis rats [1], the fungal siderophores ferricrocin (FC) and triacetylfulvarinine C (TAFC) were monitored in rat urine before and after lung inoculation with *A. fumigatus* conidia. Molecular biomarkers were separated in high/low dose infection models by high performance liquid chromatography (HPLC) and detected by mass spectrometry (MS). Rat results were compared with those in a set of 11 human clinical samples.

Results

In high dose rat model, the initial FC signal reflecting the aspergillosis appeared as early as four hours post-infection. results for seven biological replicates showed increasing exponential metabolite profiles in three-day frame of the experiment. TAFC was shown to be less populated biomarker produced by *A. fumigatus* 1059 CCF strain exhibiting identical kinetic profile to that of FC. Among nine biological replicates in the low dose model, three animals did not develop any infection up to experiment termination (day 10). One animal experienced an exponential increase of metabolites and died on the day six post-infection. All remaining five animals revealed constant or random FC levels and exhibited little or no symptoms. μ PET/CT imaging, although it was not as an early tool, further determined the infection location in vivo and allowed to visually follow the infection progression in time. In a set of 11 human ICU patients with probable aspergillosis, siderophore-based approach was more sensitive than standard galactomannan (GM) testing.

Conclusions

Siderophore detection in urine (rat or human) by mass spectrometry represents an early, innovative and

non-invasive tool for diagnosing aspergillosis. Further studies are needed to clarify the effect of antifungal prophylaxis to *Aspergillus* secondary metabolite production in human host. The general application of microbial siderophores in human diagnostics can soon represent the same clinical breakthrough like ribosomal protein typing 10 years ago.

Novel Aspect

Siderophore *Aspergillus* detection is non-invasive, rapid, sensitive and as such has a great potential to standardization in human diagnostics.

References

1. Luptakova, D., Pluhacek, T., Petrik, M., Novak, J., Palyzova, A., Sokolova, L., Lemr, K., Havlicek, V., Scientific Reports, 7, 16523 (2017).

MP-138 / ANALYSIS OF UNDERIVATIZED AMINO ACIDS IN EXHALED BREATH CONDENSATE SAMPLES AS MONITORING THERAPY FOR CHILDREN WITH TYPE 1 DIABETES

Magdalena Pyszka⁽¹⁾ - *Lucyna Konieczna*⁽¹⁾ - *Magdalena Okońska*⁽¹⁾ - *Łukasz Szymgel*⁽²⁾ - *Małgorzata Myśliwiec*⁽²⁾ - *Tomasz Bączek*⁽¹⁾

Department of Pharmaceutical Chemistry, Medical University of Gdańsk, Poland (1) - Department of Family Medicine, Medical University of Gdańsk, Poland (2)

Keywords: exhaled breath condensate (EBC) samples, amino acids, underivatization step, type 1 diabetes, LC-MS/MS

Introduction

Insulin deficiency in type 1 diabetes leads to metabolic changes. Despite technological progress, early detection of the disease is problematic and checking sugar level is burden for 66% of patients. The use of exhaled breath condensates as a non-invasive sampling method might help patients. Liquid chromatography technique coupled with mass spectrometry may help in the detection of differences in amino acid profiles in patients with type 1 diabetes and healthy ones.

Methods

A procedure for preparing EBC samples with underivatized amino acids to chromatography analysis is based on vacuum evaporation. For the determination of low molecular weight amino acids, liquid chromatography coupled to mass spectrometry (LC-MS/MS) with a hydrophilic interaction chromatography (HILIC) column was used. All amino acid profiles were determined in exhaled breath condensate samples.

Results

A sensitive and specific method was used for the simultaneous determination of the entire L-amino acid profile in exhaled breath condensate samples, a matrix taken in a completely non-invasive way, in pediatric patients with type 1 diabetes, and the results were compared with the control group. The use of vacuum evaporation as a sample preparation step has reduced the consumption of solvents, which is consistent with the principles of green chemistry. Thanks to the use of the HILIC column, it was possible to identify 23 amino acids without the necessity to derivatize them. The condensates of exhaled breath of pediatric patients with diabetes in whom the disease has been diagnosed have been analyzed.

Conclusions

The results indicate that it seems to be advisable to use the amino acid analysis in EBC, which is a method that does not burden the patient doing multiple blood glucose measurements and treated with insulin. Due to the chronic nature of diabetes, chronic complications are important for the patients' survival time and

the quality of life. The EBC can be used as a new method allowing for their early detection and can positively affect the effects of treatment of type 1 diabetes.

Novel Aspect

A new research method based on LC-MS/MS has been developed, allowing for non-invasive diagnostics of metabolic changes occurring in type 1 diabetes.

References

1. Koniczna L., Pyszka M., Okońska M., Niedźwiecki M., Bączek T., J. Chromatogr. A 1542, 72-81 (2018).
2. Nagao K., Yamakado M., Curr. Opin. Clin. Nutr. Metab. Care. 19, 328–335 (2016).
3. Karyakina A. A., Nikulina S. V., Vokhmyanina D. V., Karyakina E. E., Anaev E. K. H., Chuchalin A. G., Electrochem. Commun. 83, 81–84 (2017).
4. Kuban P., Foret F., Anal. Chim. Acta 805, 1–18 (2013).
5. Wang T. J., Larson M. G., Vasan R. S. et al., Nat. Med. 17, 448 – 453 (2011).

MP-139 / DEVELOPMENT AND VALIDATION OF A HIGH SENSITIVITY ASSAY OF ESTROGENS IN HUMAN PLASMA BY UHPLC-MS/MS WITHOUT DERIVATIZATION

Mikaël Levi⁽¹⁾ - *Hironori Kobayashi*⁽²⁾ - *Jun Watanabe*⁽¹⁾

Shimadzu Corporation, MS BU, Kyoto (1) - Shimane University Faculty of Medicine, Department of Pediatrics, Shimane (2)

Keywords: Estrogens, plasma, high-sensitivity, UHPLC-MS/MS.

Introduction

Estrogens play important roles in many biological functions as primary female sex hormones. In males, estrogens are present at lower concentrations but are also involved in several mechanisms. Modification of circulating estrogen levels can be associated with both positive and negative health effects.

In the last decade, LC/MS/MS has proven its capabilities in terms of sensitivity and specificity to become the gold standard methodology in steroid analysis. However, measuring estrogens at very low level remained challenging. Many published methods involve complex sample preparation, large sample volumes or derivatization. These processes are not compatible with a routine laboratory workflow. Here we present a UHPLC-MS/MS method to measure estrone (E1), estradiol (E2) and estriol (E3) together in the low pg/mL range.

Methods

300 µL of plasma were spiked with 30 µL the ISTD solution (deuterated E1, E2 and E3) and then loaded on a SLE+400 96-well-plate (Biotage, Sweden). After 5 min incubation, target compounds were eluted with 2x 500 µL of Ethylacetate/Hexane 75/25 v/v. After evaporation and reconstitution, extracts were analyzed using a Nexera X2 coupled to a LCMS-8060 tandem mass spectrometer (Shimadzu Corp., Japan). Using a gradient of water and methanol with NH₄F, all compounds were separated on a Raptor Biphenyl column (Restek, USA) with a cycle time of 3.5 min. Data acquisition used MRM mode with the electrospray ion source operating in negative mode.

Results

To fully optimize the detection of estrone (E1), estradiol (E2) and estriol (E3) in plasma, the method considered the impact of alternative sample preparation protocols in addition to enhancing ion sampling efficiency and MS/MS ion signals.

Several sample preparation procedures were evaluated including protein precipitation (PPT), solid phase extraction (SPE), liquid-liquid extraction (LLE) and supported liquid extraction (SLE). In this work, SLE

resulted in the highest extraction recovery with a minimal matrix interference. SLE was also considered to be a simple technique to use routinely with a relatively fast sample cycle time and created further opportunities in automation. Extraction recovery was between 96 to 97% for E1, 88 to 102% for E2 and 94 to 108% for E3. Matrix effect ranged from 0.90 to 1.00 for E1, 1.19 to 1.29 for E2 and 0.91 to 1.11 for E3.

As SLE selectively extracted E1, E2 and E3 whilst minimizing possible co-extraction of isobaric matrix interferences the peak area variance for the internal standards (ISTD) in patient samples and QC was less than 16% (%RSD was 16%, 9% and 11% for D4-E1, D3-E2 and D3-E3 respectively, n=27).

As a result, we obtained limits of quantification of 0.5 pg/mL (E1), 1 pg/mL (E2) and 5 pg/mL (E3). To our knowledge, such limits of quantification have never been reported with a simple sample preparation and no derivatization. The linear calibration range was then 0.5-500 pg/mL, 1-500 pg/mL and 5-1000 pg/mL for E1, E2 and E3 respectively. Such sensitivity was successfully used to report samples that could not be quantified before with immunoassays.

The assay also generated acceptable precision and accuracy in repeated use meeting the needs of a routine steroid hormone analysis. In validating the method certified reference serums and individual patient samples with known estrogens concentrations were measured and good correlation was obtained.

Conclusions

A method for routine high-sensitivity analysis was successfully developed for estrogen analysis in human plasma. The sample preparation is easy and very efficient to ensure low pg/mL limits of quantification as well as ruggedness and accuracy. Linearity of the assay allows to measure all kind of samples, not only the low-level ones. This method could be helpful to support clinical research.

Novel Aspect

Development of a highly sensitive routine platform for estrogens analysis in clinical pathology laboratories.

MP-140 / ICP-MS IMMUNOASSAY USING METAL-DOPED NANOPARTICLES FOR MULTIPLEX DETECTION OF HIV-1 P24 AND HBSAG

Heung Bin Lim - Chanmi Kim

Department of Chemistry, Dankook University, #119 Dandaero, Dongnam-gu, Chunan-si, Chungnam, South Korea 31116

Keywords: ICP-MS immunoassay; Nanoparticles; Determination of viruses; Multiplex detection; Clinical analysis.

Introduction

Because pathogen species can interact within the host and coinfection is quite common, it is crucial to detect viruses simultaneously in clinical analysis for accurate and fast diagnosis. So far, we have developed ICP-MS immunoassay using synthesized nanoparticles for multiplex detection and high sensitivity through particle tagging and signal amplification. And, the method has been demonstrated in a wide range of biotargets from small molecules to large cells.

Methods

In this work, we expand the ICPMS immunoassay to HIV-1 p24 and HBsAg selected among the representative viruses. For their detection, metal-doped magnetic nanoparticles and silica nanoparticles were synthesized for target extraction and tagging as probes, respectively, and sandwich-type conjugates of selected proteins representing the targets were formed through immunoreaction. Then, the quantitative determination was performed using ICPMS.

Results

For target extraction and multiplex detection, Fe₃O₄ magnetic nanoparticles were doped with Cs, and silica

nanoparticles were doped with Gd/RhBITC for HIV p24 and Eu/RhBITC for HBsAg, and then both were immobilized with the antibodies of HIV-1 p24 and HBsAg for immunoreaction, confirmed by a confocal laser fluorescence microscope and SEM. After the formation of sandwich-type conjugates, the concentration of both targets was simultaneously determined through a ratiometric detection of the doped metals using ICP-MS. The calibration curve showed a linear range of 0.78-8.0 ng/mL and 10-100 pg/mL with a detection limit of 0.61 ng/mL and 11.54 pg/mL, respectively. The multiplex detection of HIV-1 p24 and HBsAg spiked in diluted human serum produced recovery rates of 85-95%.

Conclusions

So far, we have demonstrated ICP-MS immunoassay for a multiplex detection of biomarkers, CRP (cardiovascular disease), AFP (tumor), and NSE (heart disease). In this work, we successfully expand its application to simultaneous determination of representative viruses, HIV-1 p24 and HBsAg in human serum for the first time, which implies great potential for clinical application.

Novel Aspect

In this work, we successfully demonstrated ICP-MS immunoassay using magnetic nanoparticles to determine viruses simultaneously for the first time with high sensitivity.

MP-141 / SOLID PHASE MICROEXTRACTION COUPLED TO MASS SPECTROMETRY VIA MICROFLUIDIC OPEN INTERFACE FOR POINT-OF-CARE TESTING OF TRANEXAMIC ACID IN PLASMA OF CARDIAC SURGICAL PATIENTS.

Nikita Looby⁽¹⁾ - **Marcos Tascon**⁽¹⁾ - **Vinicius Acquaro Jr**⁽²⁾ - **Nathaly Reyes-Garcés**⁽¹⁾ - **Tijana Vasiljevic**⁽¹⁾ - **Germán Augusto Gómez-Ríos**⁽¹⁾ - **Marcin Wąsowicz**⁽³⁾ - **Janusz Pawliszyn**⁽¹⁾

University of Waterloo, University of Waterloo/ Department of Chemistry, Waterloo (1) - Universidade de São Paulo, Universidade de São Paulo/ Departamento de Química, São Paulo (2) - Toronto General Hospital, Department of Anaesthesia and Pain Management, Toronto (3)

Keywords: Point-of-care, solid phase microextraction, therapeutic drug monitoring, microfluidic open interface, ambient mass spectrometry

Introduction

Tranexamic acid (TXA) is an antifibrinolytic agent routinely used in clinical settings to stop or reduce blood loss during surgeries. Attainment of TXA blood plasma concentration within the therapeutic range is difficult because it presents high interpatient variability besides, persistently high TXA concentrations have been associated with seizures. (1) Therefore, a quantitative point-of-care method based on SPME coupled to mass spectrometry was developed and validated.

Methods

TXA was monitored directly from plasma of 6 patients undergoing cardiac surgery using biocompatible SPME (BioSPME) fibers. The coating, 4 mm-long and 25 μm thickness, used HLB® particles (≤5μm) embedded in a polyacrylonitrile polymer. A microfluidic open interface (MOI) was used to couple these fibers directly to MS. (2) The developed point-of-care method was cross validated against a SPME-LC-MS/MS method using Passing Bablok regression and Bland-Altman plot.

Results

BioSPME-MOI-MS/MS revealed that the point-of-care method herein developed was able to accurately reproduce the profiles of 6 cardiac surgical patients as the profiles did not show statistical difference to those generated by the bioSPME-LC-MS/MS method. Statistical cross validation by means of the Bland-Altman plot and the Passing Bablok regression yielded for the former that 92% of the data pairs compared

were distributed randomly within the acceptable range of the limits of agreement, demonstrating that the method had no biased tendencies for the results. The latter supported that the two methods demonstrated good statistical agreement achieving a slope of 0.984. Given that the expected concentration range of tranexamic acid is between 80 – 500 µg/mL, the method produced an acceptable linear dynamic range of 50-1000 µg/mL. One of the main advantages of this approach involves a single sample turn-around time of 15 min while in a high-throughput mode (96-well plate), the speed of analysis is below 30 seconds.

Conclusions

Bio-SPME-MOI-MS/MS proves to be a reliable and feasible tool for rapid therapeutic drug monitoring. Affording total times of analysis per sample as low as 30 seconds in its high-throughput mode configuration, the currently presented method offers comparable results to the gold standard HPLC-MS/MS method with acceptable RSDs, while only necessitating a fraction of the time needed for implementation of the HPLC-MS/MS method

Novel Aspect

High throughput therapeutic drug monitoring achieved in a fraction of the time via point of care SPME-MOI-MS/MS.

References

- 1 I. Lecker, D. S. Wang, P. D. Whissell, S. Avramescu, C. D. Mazer and B. A. Orser, *Ann. Neurol.*, 2016, 79, 18–26.
- 2 M. Tascon, M. N. Alam, G. A. Gómez-Ríos and J. Pawliszyn, *Anal. Chem.*, 2018, 90, 2631–2638.

MP-142 / STUDY OF NON-COVALENT INTERACTION BETWEEN ALS-LINKED MUTANT PFN1 PROTEIN AND NATURAL COMPOUND BY MASS SPECTROMETRY

Mina Nekouei⁽¹⁾ - *Parviz Ghezellou*⁽¹⁾ - *Mahmoud Kiaei*⁽²⁾ - *Alireza Ghassempour*⁽¹⁾

Shahid Beheshti University, Medicinal Plants and Drugs Research Institute, Tehran (1) - University of Arkansas for Medical Sciences, Department of Pharmacology and Toxicology, Arkansas (2)

Keywords: Profilin 1, Protein aggregation, bottom-up proteomics, protein-ligand interaction

Introduction

Profilin 1, encoded by the PFN1 gene as one of the recently identified ALS genes, is an evolutionarily highly conserved protein[1]. Mutant PFN1 proteins form aggregate in vivo and in vitro which can be a cause of motor neuron death[2]. Mass spectrometry-based techniques can open new insights into characterization of biomarker proteins in neurodegenerative disorders. Hence, we have investigated PFN1G118V aggregation and its inhibition using mass spectrometry.

Methods

Proteins in aggregates were identified by bottom-up proteomic analysis using LC-MS/MS and MALDI-TOF. The effect of purified natural compounds was investigated by western blotting in transiently transfected (PFN1G118V and PFN1WT) cells. Moreover, we expressed PFN1WT and ALS-linked PFN1G118V in *E. coli* and purified them using anion-exchange chromatography. Recombinant PFN1G118V was incubated with compounds at various ratios and analyzed by mass spectrometry.

Results

The primary result of bottom-up proteomic analysis on the insoluble fraction of spinal cord lysates of hPFN1G118V transgenic mice revealed the identity of PFN1 along with other up- and down-regulated proteins in aggregates. Additionally, our data on the insoluble fractions of transfected and non-transfected

cells showed that only mutant PFN1 forms aggregate compared with wild type and non-transfected cells. In order to block the aggregation of mutant PFN1, we tested natural compounds derived from a medicinal plant and observed that treatment with our ROG-102 compound decreased the aggregation of PFN1G118V in vitro. To test whether this compound can directly interact with PFN1G118V and inhibit its aggregation, we sought to assess it in a cell-free system. The non-covalent interaction of recombinant PFN1G118V with the selected compound is being investigated kinetically by LC-MS/MS in non-reduced form and in digested form. Additionally, our result has been used for the binding site determination.

Conclusions

In conclusion, our mass spectrometry results can be a significant step on the road of detailed deciphering the neurodegenerative process. The identification of proteins in the aggregates will shed light on cellular protein interactions with mutant PFN1. Moreover, anti-aggregation properties of natural compounds may be highly valuable for drug discovery and therapeutic development for ALS. .

Novel Aspect

We were able to show the anti-aggregative activity of a natural compound on mutant PFN1. This novel finding has the potential as new therapy to treat motor neuron disease.

References

1. Wu C. H., Nature, 488 (7412), 499-503 (2012).
2. Fil D., Human molecular genetics, 26 (4), 686-701 (2017).

MP-143 / DETECTION AND QUANTIFICATION OF CARBOHYDRATE DEFICIENT TRANSFERRIN BY MALDI-COMPATIBLE PROTEIN CHIPS PREPARED BY AMBIENT ION SOFT LANDING

Petr Novak⁽¹⁾ - **Petra Darebna**⁽¹⁾ - **Michael Volny**⁽²⁾ - **Viktor Ruzicka**⁽³⁾ - **Jan Spicka**⁽⁴⁾ - **Radek Kucera**⁽⁵⁾ - **Petr Pompach**⁽⁶⁾

Institute of Microbiology, Laboratory of Structural Biology, Prague (1) - AffiPro, R&D, Mratin (2) - BioVendor, R&D, Brno (3) - University Hospital Kralovske Vinohrady, Department of Laboratory Diagnostics, Prague (4) - University Hospital in Pilsen, Department of Immunochemistry, Pilsen (5) - Institute of Microbiology, Laboratory of Structural Biology, Vestec (6)

Keywords: Ambient ion landing, Protein chips, Immuno-affinity enrichment, MALDI, Carbohydrate Deficient Transferrin

Introduction

Transferrin is the most important iron-transport carrier in the human body. Severe alcohol consumption leads to alterations in glycosylation of transferrin (CDT). Mass spectrometry can provide fast detection and quantification of transferrin isoforms since they have different molecular masses. In this study we use antibody chips prepared by ambient soft ion landing in combination with MALDI -TOF for detection and quantification of transferrin isoforms.

Methods

Protein chips were prepared by functionalization of indium tin oxide glass using anti-transferrin antibody. Several microliters of patient serum were applied on the antibody-modified spots and after incubation, washing and matrix deposition, the transferrin molecule forms were detected by MALDI-TOF mass spectrometer. The CDT values obtained by MALDI chip method were compared with the results obtained by capillary electrophoresis (CE).

Results

Functionalized MALDI chips were successfully used for enrichment and detection of CDT from human

serum. A cohort of samples from 186 patients was analyzed using functionalized MALDI chips and the results were compared with data obtained by CE. Out of the 186 samples, 44 were positively identified as belonging to alcoholic patients, while 142 were negative by MALDI chip approach.

Conclusions

Functionalized MALDI chips modified by anti-transferrin antibody prepared by ambient soft ion landing were successfully used for detection and quantification of CDT from human sera. The correlation of the data obtained by the CE and the chip-based MALDI was $r=0.986$, 95% confidence interval.

Novel Aspect

MALDI chips prepared by ion soft landing are suitable for quantification of CDT.

For information please contact: pnovak@biomed.cas.cz

MP-144 / ACHETYLCHOLINE AND ACHETYLCHOLINESTRASE IN A NON-NEURONAL TISSUE (CORNEAL EPITHELIUM)

Inger Oulie⁽¹⁾ - *Leon Reubsaet*⁽¹⁾ - *Trine Grønhaug Halvorsen*⁽¹⁾ - *Amund Ringvold*⁽²⁾

University of Oslo, School of Pharmacy, Dept. of Pharmaceutical Chemistry, Pharmaceutical analysis, Oslo (1) - National Hospital, University of Oslo, Eye Department, Oslo (2)

Keywords: Acetylcholine, Acetylcholinesterase, Corneal Epithelium Cells, Diurnal, Twilight and Nocturnal Mammals

Introduction

Acetylcholine (ACh) and Acetylcholinesterase (AChE) are integral parts of the central and peripheral nervous system. The substances have in addition been observed in tissues outside the nervous system such as corneal epithelial cells (CEC). ACh and AChE were analysed in CEC of diurnal, twilight and nocturnal mammals by tandem HPLC/MS-MS and fluorometric methods, respectively.

Methods

CEC from 6 different species (cattle, sheep, pig, horse, goat and cat) were obtained from local slaughterhouses. The CEC was rubbed off the eyes with a glass slide and homogenized using a pellet mixer for 2 min in mobile phase or MilliQ-water. For analysis of ACh, a HPLC/MS-MS method was used. For analysis of AChE, Amplite™ Fluorometric Acetylcholinesterase Assay Kit from AAT Bioquest® Inc was used according to the manufactures recommendations.

Results

The method for analysing ACh in CEC in this study was a modified version of the method described in [1]. Endogenous quaternary ammonium compounds (including ACh) were confirmed to be found in the CEC from all species as described in references [1,2]. The concentration of ACh and AChE in the CEC were varying between the different species. Diurnal species revealed a concentration of ACh in the range 0.18 to 2.80 µg ACh/mg CEC, in contrast to twilight and nocturnal species where ACh was close to zero, similar to previous observations [3]. Twilight and nocturnal species revealed a high concentration of AChE in the CEC, in contrast to diurnal species that revealed a low concentration. For AChE the concentration was varying from 0.5 mU AChE/mg CEC in cattle and horse to 20.4 mU AChE/mg CEC in cat.

Conclusions

Diurnal animals (cattle, horse goat and sheep) revealed high amounts of ACh in the CEC, in contrast to cell samples from twilight and nocturnal animals (pig and cat) that contain low amounts of ACh in CEC. For

AChE, the picture is reversed: diurnal animals revealed low amounts of AChE and twilight and nocturnal revealed high amounts of AChE in CEC.

Novel Aspect

To our knowledge, this is the first time Acetylcholine and Acetylcholinesterase has been measured in CEC from different mammals.

References

1. Reubsaet J. L. E., Ahlsen, E., Haneborg, K. G., Ringvold, A., Journal of Chromatographic Sciences, Vol. 41, 151-156 (2003).
2. Reubsaet J. L. E., Ringvold, A., Journal of Chromatographic Sciences, Vol. 43, 401-405 (2005).
3. Ringvold A., Reubsaet J. L. E., Cornea, Vol. 24, 1000-1003 (2005).

MP-145 / DEVELOPMENT AND VALIDATION OF A SENSITIVE LC-LINEAR ION TRAP METHOD FOR THE DETERMINATION OF TRACE DI-IODOTHYRONINES IN HUMAN BLOOD SERUM

Keith Richards⁽¹⁾ - ***Josef Köhrle***⁽¹⁾ - ***Daniel Rathmann***⁽²⁾ - ***Raymond Monk***⁽¹⁾

Charite Universitätsmedizin Berlin, Institute for Experimental Endocrinology, Berlin (1) - Chromicent GmbH, Pharmaceutical Development, Berlin (2)

Keywords: Thyroid hormone metabolites, serum, tandem mass spectrometry, MS3, validation

Introduction

In addition to the thyroid hormones (TH) 3,5,3'-triiodothyronine (T3) and L-thyroxine (T4), thyroid hormone metabolites (THM) exert their own distinct biological actions [1]. By incorporating MRM/MS3 switched scanning and extract de-lipidation, we have modified and validated our tandem mass spectrometry method [2] for 3,5-diiodothyronine (3,5-T2) and 3,3'-diiodothyronine (3,3'-T2) in human blood serum, and applied it to measure T2 concentrations in healthy individuals.

Methods

The method includes mixed mode cation exchange extraction (SPE) of 150µL serum, followed by de-lipidation using Waters Ostro sorbent (both using a 96-well plate format). The liquid chromatography-mass spectrometric method uses switched mode (MRM/MS3) acquisition using a linear ion trap (LIT).

Results

Endogenous 3,5-T2 and 3,3'-T2 were detected in commercially obtained serum, in addition to T3,3,5',3'-triiodothyronine (rT3) and T4. Inter- and intra-day accuracy and precision for exogenously spiked TH/THM were in the range 70-100% and less than 20% RSD respectively for low, medium and high concentrations. Excellent linearity was obtained across a range from 0.01 nM to 250 nM. Sera from healthy individuals contained concentrations of T3 and T4 within expected reference ranges (1.3 - 3.1 and 65.6-181.5 nM respectively) [3] whilst concentrations of 3,5-T2 and 3,3'-T2 were significantly lower (0.02 – 0.05 and 0.1 – 0.25 nM respectively). Limits of quantification (LOQ) for T2 isomers were 0.01 nM. Matrix effects for all analytes were 30% or less.

Conclusions

The use of switched MRM/MS3 scanning enables quantification of endogenous TH/THM with widely differing concentrations in a single analysis, with excellent inter- and intra-day accuracy and precision. Due to enhanced selectivity and low chemical noise, MS3 scanning delivered 4-5 times lower LOQ compared with traditional MRM for T2 isomers. In healthy individuals both 3,5-T2 and 3,3'-T2 appear to be minor metabolites (derived from T3 and T3/rT3 respectively).

Novel Aspect

The direct SPE extraction of serum without preliminary protein precipitation, combined with solid phase de-lipidation and switched MRM/MS3 scanning.

References

1. Jonas W., Lietzow J., Wohlgemuth F., Hoefig C. S., Wiedmer P., Schweizer U., Köhrle J., Schurmann A., Endocrinology, 156,389-399 (2015).
2. Richards K. H., Rathmann D., Monk R., Köhrle J., European Thyroid Journal, 6 (suppl 1), 57 (2017).
3. Mutlu M. et al. J Matern Fetal Neonatal Med., 25(2), 120-4 (2012).

MP-146 / USING DATA DEPENDENT ACQUISITION TO PRODUCE SPECIFIC AND RELIABLE MRM TRANSITIONS FOR SUSPECTS ON A QTRAP INSTRUMENT

Stine Krogh Steffensen - Inge Fomsgaard

Department of Agroecology, Aarhus University, Aarhus University/Department of Agroecology, Flakkebjerg

Keywords: Data dependent acquisition, HPLC- QTrap MS, Benzoxazinoids, Suspects, semi-quantification

Introduction

In complex biological samples similar small molecules often exists as different structural isomers or similar phase II conjugates such as glucosides or glucuronides. This often results in multi reaction monitoring (MRM) transitions that are too general for reliable assignment of selected unknowns or suspects in a triple quadrupole instrument. However, by letting data dependent acquisition techniques guide selection of MRMs, a reliable and sensitive method was developed.

Methods

Benzoxazinoid-rich complex biological samples were fractionated using consecutive SPE and HPLC to concentrate benzoxazinoid metabolites. Our established HPLC-MSMS (ABSciex Qtrap 4500) benzoxazinoid quantification method was adapted for data dependent acquisition (DDA/IDA) using mainly enhanced precursor scan (EMS) and MRM in different forms as survey scans and enhanced product ion (EPI) scans as dependent scans.

Results

24 different peaks were deduced as interesting based on their spectral properties and varying similarities to known benzoxazinoids or metabolites thereof. Figures will show selected examples of these such as the spectrum the glucuronide of the lactam HBOA, known from rye. This spectrum contains the fingerprint of HBOA (164/136/118/108) and of the glucuronic acid (193/175/113) as well at the neutral loss of 18 u from the precursor ion (340) typical for carbohydrates. 340/164 and 340/108 were reliable MRM transitions for this ion, however general in their nature, and aided the discovery of other interesting peaks at other retention times. 340/118 and 340/136 were also reliable but much more selective. Other examples are the MRMs 354/164 and 354/108 that gave a peak at Rt=10.9 min. However, the spectrum varied slightly from sample to sample. Based on spectral data obtained by DDA, the MRMs 354/165 and 354/118 were introduced, revealing the presence of two distinct peaks differing by only 0.05min in retention time.

Conclusions

Highly selective MRM transitions were produced that enabled reliable detection and assignment of selected unknowns and suspects in complex biological matrices. Using this DDA-guided method for selection of MRM transitions, it was possible to distinguish even between compounds with similar retention times, precursor ions, and even fragment ions.

Novel Aspect

Using DDA techniques for selection and evaluation of MRM transitions for semi-quantitative methods for selected unknowns and suspects.

MP-147 / EVALUATION OF UDP-GLUCURONOSYLTRANSFERASE 1A1 ACTIVITY TOWARDS BILIRUBIN**Atsuko Takeuchi**

Instrumental Analysis Center, Kobe Pharmaceutical University, Kobe

Keywords: UDP-glucuronosyltransferase 1A1, UGT1A1, LC/MS/MS, bilirubin, enzyme activity

Introduction

UDP-glucuronosyltransferase 1A1 (UGT1A1) catalyzes glucuronidation of bilirubin. Mutations in the UGT1A1 gene cause bilirubin glucuronidation deficiencies, leading to diseases with hyperbilirubinemia or jaundice, Crigler-Najjar syndrome type I, type II and Gilbert syndrome. Since all the patients carrying these mutations showed hyperbilirubinemia, the enzyme activities of the mutated enzymes were determined.

Methods

We applied the liquid chromatography-mass spectrometry (LC-MS) method to evaluation of UGT1A1 activity toward bilirubin.

Chromatographic separations were performed using an Ultimate 3000 HPLC system (Dionex). Mass spectrometry was performed on an LTQ Orbitrap Discovery mass spectrometer (Thermo Scientific). The electrospray ionization source was used.

Results

We determined the remaining UGT1A1 activity of the models of wild-type and mutant-type UGT1A1. The UGT1A1 activities of the models were considered as the normal and diseases with hyperbilirubinemia UGT1A1 activity, respectively.

We applied LC/MS to evaluation of mutation effects on the UGT1A1 activities towards bilirubin. Bilirubin has two sites capable of glucuronidation and forms bilirubin mono- and di-glucuronides. LC/MS can detect a small amount of their glucuronides with accuracy. We confirmed the reduced substrate-glucuronidation activity of mutated UGT1A1 enzymes.

Conclusions

We aimed to establish a highly sensitive method to determine the UGT1A1 activities. The activities varied according to wild-type and mutant-type UGT1A1. We tested here all the mutations found in patients with hyperbilirubinemia. The obtained results allow us to discuss the relationship between remaining activity and clinical phenotype of the patients with hyperbilirubinemia.

Novel Aspect

We established a highly sensitive method to determine the UGT1A1 activities and evaluated wild-type and mutant-type UGT1A1 activities.

MP-148 / A BIOMPHALARIA GLABRATA PEPTIDE THAT STIMULATES EXTREME BEHAVIOUR MODIFICATIONS IN AQUATIC FREE-LIVING SCHISTOSOMA MANSONI MIRACIDIA

Tianfang Wang⁽¹⁾ - **Russell Wyeth**⁽²⁾ - **Di Liang**⁽¹⁾ - **Utpal Bose**⁽¹⁾ - **Guoying Ni**⁽³⁾ - **Donald McManus**⁽⁴⁾ - **Scott Cummins**⁽¹⁾

University of the Sunshine Coast, School of Science and Engineering, Sippy Downs (1) - St. Francis Xavier

University, Department of Biology, Antigonish (2) - Griffith University, School of Medical Science, Gold Coast (3) - University of Queensland, QIMR Berghofer Medical Research Institute Berghofer, Brisbane (4)

Keywords: *Schistosoma mansoni*, *Biomphalaria glabrata*, miracidia, kairomone, proteomics

Introduction

The human disease schistosomiasis is caused by the parasite *Schistosoma mansoni*, which requires an intermediate host, the freshwater snail *Biomphalaria glabrata*. An excellent olfactory sense allows the parasite to locate an appropriate host. Towards identifying the kairomone biomolecule(s) that helps *S. mansoni* miracidia detect its host, we have performed an integrated analysis of the snail-conditioned water (SCW), through chromatography and bioassays of behaviour, followed by mass spectrometry.

Methods

The lyophilised SCW crude extract was resuspended in 0.1% TFA and fractionated using HPLC, and the fractions were collected and lyophilised for bioassays, with behaviour analysed by TrackMate plugin (1) and MtrackJ (2) for FIJI software (3). The LC-MS/MS analysis of metabolites and peptides were carried out using UHPLC-QToF and Nano HPLC-ESI-Triple TOF, respectively. Data analysis of was performed using Agilent MassHunter Qualitative software and PEAKs studio. CD experiments were performed in a Chirascan CD spectrophotometer (Applied Photophysics).

Results

Using an in vivo bioassay in conjunction with behaviour analysis, we first determined a raw mucus extract and filtered SCW could stimulate behaviour changes in *S. mansoni* miracidia, which was in alignment with previously classified miracidial responses (4,5). The lyophilised crude extracts from SCW were HPLC-fractionated and collected for subsequent bioassays, and a single HPLC fraction demonstrated behaviour modification bioactivity. Metabolites and proteins present within the active HPLC fraction and two adjacent fractions were identified using LC-MS/MS. We found 16 metabolites, all of which did not match to any known compound within the METLIN (6) or HMDB (7) metabolite databases. Also in the active HPLC fraction, a total of 24 peptides were identified that matched to precursor proteins in the genome-derived *B. glabrata* protein database (8). The activity analysis of these 24 peptides identified a single peptide with alpha helical characteristics stimulates extreme behaviour modifications in miracidia, and it can be produced within numerous snail tissues.

Conclusions

We have reported the identification of a peptide, which is secreted from adult *B. glabrata* and triggers extreme behaviour modifications in the *S. mansoni* miracidia. Its precursor appears to be unique to this snail, thus providing an ideal species-specific attraction cue for *S. mansoni* miracidia. This finding contributes greatly to our understanding of a key part within the parasite's life-cycle, and may help towards establishing novel control interventions for schistosomiasis.

Novel Aspect

This information will be helpful for developing approaches to manipulate this parasite's life cycle, and opens up new avenues for exploring other parasitic diseases.

References

1. Tinevez, J. Y. et al. TrackMate: An open and extensible platform for single-particle tracking. *methods* (San Diego, Calif.)115, 80-90, doi:10.1016/j.ymeth.2016.09.016 (2017).
2. Meijering, E., Dzyubachyk, O. & Smal, I. methods for cell and particle tracking. *methods in enzymology*504, 183-200, doi:10.1016/b978-0-12-391857-4.00009-4 (2012).
3. Schindelin, J. et al. Fiji: an open-source platform for biological-image analysis. *Nature Methods* 9, 676-682, doi:10.1038/nmeth.2019 (2012).

4. MacInnis, A. J. Responses of *Schistosoma mansoni* miracidia to chemical attractants. *The Journal of parasitology* 51, 731-746 (1965).
5. Haas, W., Gui, M., Haberl, B. & Strobel, M. Miracidia of *Schistosoma japonicum*: approach and attachment to the snail host. *The Journal of parasitology* 77, 509-513 (1991).
6. Smith, C. A. et al. METLIN: a metabolite mass spectral database. *Therapeutic drug monitoring* 27, 747-751 (2005).
7. Wishart, D. S. et al. HMDB 3.0-The Human Metabolome Database in 2013. *Nucleic acids research* 41, D801-807, doi:10.1093/nar/gks1065 (2013).
8. Adema, C. M. et al. Whole genome analysis of *Biomphalaria glabrata* (Mollusca, Lophotrochozoa), a snail intermediate host for transmission of schistosomiasis. *Nature communications*, in press (2016).
9. Mitta, G., Adema, C. M., Gourbal, B., Loker, E. S. & Theron, A. Compatibility polymorphism in snail/schistosome interactions: From field to theory to molecular mechanisms. *Developmental and comparative immunology* 37, 1-8, doi:10.1016/j.dci.2011.09.002 (2012).

MP-149 / ANALYSIS OF SWEETENERS IN URINE USING NOVEL STATIONARY PHASE AND SELECTIVE LC-MS/MS METHOD

Barbora Cervinkova⁽¹⁾ - ***Richard Kay***⁽²⁾ - ***Albert Koulman***⁽³⁾

NIHR BRC Nutritional Biomarker Laboratory, University of Cambridge, Cambridge, United Kingdom (1) - Metabolic Research Laboratories, University of Cambridge, Cambridge, United Kingdom (2) - NIHR BRC Nutritional Biomarker Laboratory, Metabolic Research Laboratories, University of Cambridge, Cambridge, United Kingdom (3)

Keywords: sweeteners, urine, biomarkers, HILIC, UPLC-MS/MS

Introduction

There are only a limited number of published methods for monitoring of sweeteners intake. Most of these are focused on determination of small biomarker groups [1] and investigation their role in human metabolism. Moreover, these methods were either published deep in the past[2,3] or use chemical derivatisation [4,5]. In our study we aim to develop a novel and simple LC-MS/MS method for measurement of several intake biomarkers for different groups of sweeteners.

Methods

Analysis of the sweeteners was performed using an ACQUITY UPLC System (Waters, Elstree, United Kingdom) coupled to a Xevo TQ-XS Triple Quadrupole mass spectrometer (Waters, Elstree, United Kingdom). For optimal separation of structurally related compounds different stationary phases, including C18, Amide, HILIC and Biphenyl, were tested. Moreover, in order to achieve greater ionization efficiency for small molecules novel UniSpray ion source was tested.

Results

Final separation was achieved using a HILIC stationary phase with gradient of buffer and acetonitrile. The mobile phase gradient separated artificial sweeteners from sugar alcohols and natural sugars over a total 7.0 minute run time. For accurate detection of analytes ionization in negative mode and specific MRM transitions were used.

Conclusions

The novel UPLC-MS/MS method with HILIC separation of structurally related sweeteners in human urine was developed in this study. This method was sensitive enough to comprehensively detect the broad spectrum of artificial and natural sweeteners without applying complicated derivatisation steps.

Novel Aspect

This unique method enables rapid analysis of several different groups of sweeteners in human urine in one run. It may serve to identify biomarkers in nutritional research.

References

1. Logue C., Dowe L. R., Strain J. J., Verhagen H., McClean S., Gallagher A. M., *Journal of Agricultural and Food Chemistry*, 65, 4516-4525 (2017).
2. Wilson S. A., Wilkinson K., Crews H. M., Davies A. M., Dick C. S., Dumsday V. L., *Food Additives & Contaminants*, 16, 227-238 (1999).
3. Bekersky D., Colburn W. A., *Analytical Letters*, 13, 805-816 (1980).
4. Licea-Perez H., Jonnotula V., Zohrabian S., Karlinsey M., *Analytical Methods*, 8, 3023-3033 (2016).
5. Casals I., Reixach M., Amat J., Fuentes M., Serra-Majem L., *Journal of Chromatography A*, 750, 397-402 (1996).

MP-150 / PHARMACOKINETICS AND BIOAVAILABILITY OF DIFFERENT ACETYLSALICYLIC ACID FORMULATIONS ASSESSED BY LIQUID CHROMATOGRAPHY-TANDEM MASS SPECTROMETRY IN HEALTHY SUBJECTS

Jessica Rizzo - Michele Dei Cas - Federico Maria Rubino - Rita Clara Paroni

University of Milan, Department of Health Science, Milan, Italia

Keywords: pharmacokinetics, LC-MS/MS, aspirin, enteric-coated formulation

Introduction

Low-dose acetylsalicylic acid (ASA, 100 mg/die) is used in thromboprophylaxis. Enteric-coated formulation (EC-ASA) is commonly used for its lower risk of side effects. Some patients on EC-ASA do not respond appropriately and recent studies showed that poor responsiveness is more frequent with EC-ASA [1]. Aim of this study was to validate a method useful to study the pharmacokinetics (PK) of ASA in healthy subjects treated with two different aspirin formulations.

Methods

Liquid chromatography tandem mass spectrometry (LC-MS/MS) technique was used for determination of ASA and salicylic acid (SA) in human plasma, using the respective deuterated isotopomers as internal standards (ASA-d₄, SA-d₄). The method was validated according to FDA guidelines and was applied to evaluate the PK of ASA and SA at different times (30 min-24 h) in healthy subjects (n=11) treated for 1 week with 100 mg of EC-ASA or plain Aspirin (plain-ASA).

Results

Analytes and internal standards were extracted from human plasma by protein. The compounds were separated on a reverse-phase column with an isocratic mobile phase. Compounds parameter were optimized to have the higher signal to noise ratio in negative multiple reaction mode. Ion transition recorded were 179→137 for ASA, 183→141 for ASA-d₄, 137→93 for SA and 137→141 for SA-d₄. The calibration curves were linear ($r^2 \geq 0.99$) over the concentration range of 20-2000 ng/mL for ASA and 20-8000 ng/mL for SA. LOQ was 20 ng/mL for both.

Plasma AU Cafter EC-ASA administration was 702(604-881) ng h/mL (median, 25%-75% CI). t_{max} occurred between 2-8 hours after intake. Two subjects did not show ASA absorption within the observation period confirmed by SA trend. PK of plain-ASA showed AUC of 823 (635-1013) ng h/mL and all subjects absorbed the drug with t_{max} between 0.5-1 h.

Conclusions

The developed method allowed to study the analytical fate of aspirin "in vivo" by comparing two different

ASA formulations. Coated formulation showed a variable behavior and in 2 subjects absorption seemed impaired. On the contrary using plain-ASA variation is reduced and all subjects absorbed efficiently the drug. Causes of inadequate response to ASA may be related to different gastro-intestinal availability of coated formulation.

Novel Aspect

Availability of a validated method is pivotal to study ASA absorption and metabolism in different cohorts of patients known to have poor pharmacodynamics responses to the drug.

References

1. Grosser T, Fries S, Lawson JA, Kapoor SC, Grant GR, FitzGerald GA, Circulation , 127, 377-85 (2013).

MP-151 / WHICH LC-MS/MS PLATFORM IS MOST APPROPRIATE FOR THE QUANTITATIVE ANALYSIS OF STEROIDS IN URINE, SERUM AND ORAL FLUID FOR CLINICAL RESEARCH

Sergio Indelicato⁽¹⁾ - *Rory M Doyle*⁽²⁾ - *Douglas McDowell*⁽²⁾

Thermo Fisher, CMD, Paris, France (1) - Thermo Fisher, CMD, Somerset, New Jersey, United States (2)

Keywords: steroids, serum, triple quadrupole, liquid-liquid extraction

Introduction

Steroid hormones control many physiological aspects in our body. A robust, sensitive, accurate and specific LC/MS/MS analytical method has been developed for the quantitation of 19 steroids in serum. A simple liquid-liquid extraction sample preparation allows achieving the required sensitivity. This research study was carried out to evaluate which LC-MS/MS instrument is best suited to characterize and quantitate the steroids in serum.

Methods

Thermo Scientific™ TSQ Fortis™, Thermo Scientific™ TSQ Quantis™, and Thermo Scientific™ TSQ Altis™ tandem mass spectrometers and a Thermo Scientific™ Vanquish™ HPLC system were utilized. A 200 µL serum sample was used and a Thermo Scientific™ Accucore™ C18 column achieved baseline chromatographic separation. Quantitative analysis was performed using SRM scan in positive and negative mode. Accuracy of the analytical method was verified using reference samples.

Results

Good linearity and reproducibility were obtained with the concentration range from 10 pg/ml to 1000 ng/ml for the steroids with a coefficient of determination >0.995 for all mass spectrometer platforms. The lower limit of detection (LLOD) and lower limit of quantitation (LLOQ) were determined. Excellent reproducibility was observed for all platforms, being the %CV for each calibration point always < 10% for all steroids. The inter-assay precision was determined by extracting and quantifying three replicates of in-house control material resulting in %CV for the steroids of < 10% deviation from the targeted mean. The inter-assay precision was determined over three consecutive days and was found to have a %CV < 10% for each steroid within its respective linear range for the three levels of pooled serum sample control material respectively. The accuracy was determined by the analysis of in-house control material as the percentage deviation from the targeted mean. The results were < 10% for all levels in each matrix.

Conclusions

A sensitive, simple, specific and accurate liquid chromatography QQQ mass spectrometry analytical method was developed and verified for the measurement of 19 steroids in serum and applied to any Thermo Scientific™ tandem mass spectrometer platform. The sample preparation technique is quick and easily

applicable for high throughput analysis. The reported method achieves research laboratory required analytical performances for steroids analysis in serum.

Novel Aspect

Utilization of various triple quadrupole mass spectrometers to determine the most appropriate platform for steroid analysis in serum.

MP-152 / GC-MS AND LC-MS PROFILING OF EXHALED BREATH CONDENSATE: INVESTIGATION OF SAMPLING CONSIDERATIONS

Anna Karen Carrasco Laserna⁽¹⁾ - *Hui Fang Lim*⁽²⁾ - *Nadia Suray Tan*⁽²⁾ - *Elvi Horiyanto*⁽¹⁾ - *Bao Hui Ng*⁽¹⁾ - *Nicholas Paton*⁽³⁾ - *Sam F. Y. Li*⁽¹⁾

National University of Singapore, Department of Chemistry, Singapore (1) - National University Hospital, Respiratory and Critical Care Medicine, Singapore (2) - National University Hospital, Infectious Diseases Department, Singapore (3)

Keywords: exhaled breath, exhaled breath condensate, metabolomics, VOCs, sputum induction

Introduction

Exhaled breath condensate (EBC) is a promising specimen for use in diagnostics and therapy monitoring. It consists of condensed water vapor, water soluble volatile organic compounds (VOCs), as well as non-volatile hydrophilic and hydrophobic components [1]. It can be easily and non-invasively obtained, presenting an advantage over sputum, which is a more commonly used sample for respiratory disease assays, but is hard to obtain from children and weak patients.

Methods

In this study, we explore the effects of meal and sputum induction (SI) on the EBC profile. Samples were collected from 10 healthy volunteers, divided into 2 batches. For each batch, samples were collected for 3 days at different timepoints, before and after Sland meals. Double filtration sterilization was done in consideration of plans for infectious disease biomarkers study. VOCs were analyzed using GC-MS and the non-volatile components using LC-MS.

Results

Preliminary analysis of GC-MS and LC-MS results using principal components analysis (PCA) showed inter-batch effects for the GC-MS profiles across three days of sampling indicating possible differences in the VOC composition in the sampling environment of the two batches. It is also observed that for Batch 1, there is some separate clustering of samples before and after SI. This could be indicative of a change in the EBC VOC profile after SI. Though there were still some inter-day batch effects observed, samples before a low-fat breakfast clustered closely while the samples for each batch collected 30 minutes and 2 hours after the meal clustered together. Samples before and after a high-fat breakfast were found to cluster together. PCA of LC-MS data showed distinct clustering of samples before and after SI while samples before and after meals were mostly overlapped.

Conclusions

Results revealed that there can be changes in the VOCs and non-volatile components of exhaled breath after sputum induction. The VOCs profile can also be affected by the meal taken before sampling while the profile of non-volatile components does not significantly vary whether sampling is done before or after meal. This indicates that EBC collection should not be performed after sputum induction of patients and at least 2 hours after a meal in breath biomarkers studies.

Novel Aspect

This study demonstrates multi-platform profiling of EBC components using GC-MS and LC-MS to establish the effects of meals and sputum induction in EBC collection.

References

1. HorváthI., Hunt J., Barnes PJ., Eur Respir J, 26, 523-548 (2005).

MP-153 / SET UP OF A SIMPLE, FAST AND FLEXIBLE LC-MS/MS ID METHOD FOR THE SIMULTANEOUS DETERMINATION OF STEROIDS IN DIFFERENT BIOLOGICAL SAMPLES

Gianluca Bartolucci⁽¹⁾ - *Marco Pallecchi*⁽²⁾ - *Sebastiano Gozzo*⁽²⁾ - *Giovanna Danza*⁽²⁾

University of Florence, Department of Neuroscience, Psychology, Drug Area and Child health, Firenze, Italy (1) - University of Florence, Department of Experimental Biomedical and Clinical Science, Firenze, Italy (2)

Keywords: steroids, mass spectrometry

Introduction

LC-MS/MS ID is nowadays recognized as the gold standard method to quantify steroids and clinical laboratories are adopting MS for steroids analysis to overcome inaccuracy of immunometric methods. [1] Costs of LC-MS/MS are now competitive, but this technique is still less automatized, and methods are less standardized. [2] In setting up this new method we focused our attention on characteristics that facilitate routine application also with different instruments.

Methods

We set up a two dimensional chromatography. The first column (C18 Luna 5 μ m 2 cm x 2 mm) was used as sampling column and flushed with a polar eluent to eliminate polar compounds and concentrate steroids; the analytical column (C18 luna 3 μ m 5 cm x 2 mm) separated androgens, progestogens, glucocorticoids that were quantified with a SCIEX 6500QTRAP in positive mode. Different eluents [3] and a PFP 3 μ m 5 cm x 2 mm analytical column allowed quantification of estrogens in negative mode.

Results

This method allows the simultaneous determination of testosterone, androstendione, progesterone, 17OH progesterone, dihydrotestosterone, cortisol and cortisone, or estrone and estradiol, in a 12 min run time with high sensibility. Samples are simply diluted with the internal standard mix in methanol and water, centrifuged to eliminate proteins, and injected. It was validated calculating precision, accuracy and sensibility for each compound. To assess the method flexibility we tested it with different biological samples (urine, serum, tissue homogenates, follicular fluids), and established that the minimum sample dilution to avoid matrix effect and preserve chromatographic performance is 1:5 with an injection volume range from 10 to 100 μ L. Using the SCIEX 6500QTRAP we dilute urine and serum samples 1:10 with an injection volume of 25-50 μ L. The method is already used in our hospital (AOU-Careggi- Firenze)) for routinely analysis of urinary free cortisol and of cortisol/cortisone ratio and It will be soon adopted for serum steroids.

Conclusions

The method fit the routine characteristics being simple, fast and robust. The sensitivity reached allows an accurate determination of androgens also in women and pre-puberal subjects. The efficiency of the proposed chromatography confers flexibility necessary to the scope. In fact, the possibility to change sample dilution or injection volume, without affecting performance, are useful to adapt the method to different sample and different instruments.

Novel Aspect

The proposed LC-MS/MS system allows the analysis of steroids in many biological matrices with minimal sample preparation, maintaining high level of detection performances.

References

1. Adaway JE, Keevil BG, Owen LJ, Ann Clin Biochem. 52, 18-38, (2014)
2. Zendjabil M, Chellouai Z, Abbou O, AnnEndocrinol (Paris), 77, 43-8, (2016)
3. Fiers T, Casetta B, Bernaert B, et al., J Chromatogr B Analyt Technol Biomed Life Sci, 15, 893-894:57-62, 2012)

MP-154 / CHEMOMETRIC OPTIMIZATION OF UPLC-MS/MS ASSAY FOR CLINICAL DIAGNOSIS AND PHARMACOTHERAPY MONITORING OF PATIENTS WITH APRT DEFICIENCY

Margrét Thorsteinsdóttir⁽¹⁾ - *Unnur Arna Thorsteinsdóttir*⁽¹⁾ - *Finnur Freyr Eiriksson*⁽²⁾ - *Hrafnhildur L. Runolfsson*⁽³⁾ - *Vidar O Edvardsson*⁽⁴⁾ - *Runolfur Palsson*⁽⁴⁾

University of Iceland, Faculty of Pharmaceutical Sciences, Reykjavik, Iceland (1) - ArcticMass, Mass Facility, Reykjavik, Iceland (2) - University of Iceland, Faculty of Medicine, Reykjavik, Iceland (3) - Landspítali – The National University Hospital of Iceland, Faculty of Medicine, Reykjavik, Iceland (4)

Keywords: Clinical mass spectrometry, UPLC-MS/MS, Chemometrics, Design of experiments, APRT deficiency

Introduction

Patients with APRT deficiency have excessive renal excretion of 2,8-dihydroxyadenine (DHA), causing chronic kidney disease. Treatment with allopurinol or febuxostat reduces DHA excretion and kidney stones formation. A chemometric approach was used for optimization of a UPLC-MS/MS assay for simultaneous quantification of DHA, adenine biomarkers and the pharmacological agents in human urine and plasma for clinical diagnosis and pharmacotherapy monitoring.

Methods

Design of experiments (DoE) was used for optimization of UPLC-MS/MS assay for quantification of DHA, adenine, adenosine, inosine, hypoxanthine, xanthine and the drugs allopurinol, its active metabolite, oxypurinol, as well as febuxostat in human plasma and urine. Different sample preparation methods were tested prior to analysis of samples from controls and patients with APRT deficiency, before and after treatment with allopurinol or febuxostat.

Results

A UPLC-MS/MS assay for simultaneous quantification of DHA and the main purine metabolites in human plasma and urine was developed by chemometric approach. D-optimal design was used to reveal both quantitative and qualitative significant factors influencing retention time and peak area for the analytes. Significant factors were optimized with CCF-design and related to sensitivity utilizing PLS regression. There was a strong interaction effect between several variables, indicating that these variables cannot be independently controlled to reach optimum settings. The assay was used to measure both plasma and urine concentration of the analytes in samples from healthy controls and APRT patients, before and after treatment. The DHA concentrations in plasma from 3 APRT patients were 456, 459 and 741 ng/mL off drug treatment and 130, 61 and 27 ng/mL while on allopurinol or febuxostat therapy. DHA was not detected in plasma from healthy controls. Highly significant changes were observed in the urinary excretion of DHA and adenine with drug therapy.

Conclusions

Optimizing UPLC-MS/MS assay by DoE, only a fraction of experiments was used compared to what would

have been required by changing one separate factor at a time (COST) approach. Our data suggest that the UPLC-MS/MS assay will greatly facilitate clinical diagnosis of patients with the rare hereditary purine metabolism disorder, APRT deficiency. At the same time the clinical assay can be used for monitoring pharmacotherapy and treatment adherence among these patients.

Novel Aspect

Implementation of UPLC-MS/MS assay to facilitate clinical diagnosis and therapeutic drug monitoring of patients with the rare kidney disorder APRT-deficiency.

MPS-S07 - Quantitative proteomics

Chairs: Michael Linscheid, Nathalie Norais

MP-155 / OPTIMIZATION OF BIOTIN ENRICHMENT PROTOCOL AND ITS APPLICATION FOR THE ANALYSIS OF ORGANELLE-ENRICHED PROTEINS BY LC-MS/MS

Shiva Ahmadi - Dominic Winter

Institute of Biochemistry and Molecular Biology, University of Bonn, Bonn

Keywords: BioID, Affinity Enrichment, Organelle, MS

Introduction

The strong interaction of biotin and avidin is widely used for the enrichment of proteins from complex mixtures. The recovery of biotinylated proteins and the discrimination between unspecific binding and specifically enriched proteins remain major challenges. We compared different strategies for the purification of biotin-enriched proteins followed by MS analysis and applied them to the analysis of organelle proteomes.

Methods

Equal amounts of heavy control and light stably BioID [1] transfected SILAC cells were combined, lysed and processed with 4 methods (three replicates): 1,2: Binding to beads followed by digestion and off-gel fractionation using 2 sets of washing buffers [2]; 3. Digestion prior to biotin capture [3]; and 4. Elution from the beads and in-gel digestion [4]. Organelle proteomes were enriched using different BioID constructs. Samples were analyzed using LC-MS/MS.

Results

Using unbiotinylated 3plex SILAC labeled HEK cells (in two replicates), a significance threshold for biotinylation (95% of proteins) was calculated. Peptides with experiment/control ratios above this threshold were considered to be specifically enriched and therefore biotinylated. On bead digestion followed by offgel fractionation resulted in identification of 3165 and 3006 proteins on average for 2 different sets of washing buffers (82.3% and 81.6% biotinylated). Incubation of peptides with beads yielded on average in 173 proteins identifications (64% biotinylated) while elution of biotinylated proteins followed by in gel digestion resulted on average in 711 proteins (81.3% biotinylation). Using method 1, cells transiently transfected with 4 different BioID constructs (three replicates) located at the cytoplasm, endoplasmic reticulum, nucleus, and mitochondria were analyzed, resulting in distinct subcellular proteomes.

Conclusions

Except digestion of proteins prior to on bead biotin capture, all approaches resulted in similar enrichment

efficiencies for biotinylated proteins. On bead digestion resulted in the highest number of identifications being in our hands the best method for the analysis of affinity captured biotinylated proteins. It was possible to identify subcellular proteomes using different BioID constructs.

Novel Aspect

Comparison of methods for avidin affinity capture of biotinylated proteins for MS analysis and analysis organelle proteomes using BioID constructs.

References

1. Roux K. J., Kim D. I., Raida M., Burke B., J. Cell Biol., 196, 801–810 (2012).
2. Fukuyama H., Ndiaye S., Hoffmann J., Rossier J., Liuu S., Vinh J., Verdier Y., J. Proteomics, 75, 4610–4619 (2012).
3. Schiapparelli L. M., McClatchy D. B., Liu H. H., Sharma P., Yates J. R., Cline H. T., J. Proteome Res., 13, 3966–3978 (2014).
4. Rösli C., Rybak J., Neri D., Elia G., Methods Mol. Biol., 418, 89-100 (2008).

MP-156 / SAMPLE MULTIPLEXING BY ISOELECTRIC FOCUSING (PI-CODING) FOR ISOBARIC LABELING AND LABEL-FREE QUANTIFICATION

Juan Astorga-Wells⁽¹⁾ - **Christian Beusch**⁽¹⁾ - **Thorleif Lavold**⁽²⁾ - **Roman A. Zubarev**⁽¹⁾

Karolinska Institutet, Medical Biochemistry and Biophysics, Stockholm (1) - Biomotif AB, Stockholm (2)

Keywords: Multiplexing, Isoelectric Focusing, Isobaric Labeling

Introduction

Sample multiplexing is an important tool for improving the throughput of MS analysis. Herein, we introduce the concept of Isoelectric Point Coding (pI-Coding) for an additional x2 or x3 level-multiplexing power of isobaric labeled samples. An in-solution isoelectric focusing device for fast pI-Coding is described, especially suitable for combining biological-replicates into a single sample. The utilization of pI-Coding for LFQ and TMT samples is demonstrated.

Methods

Steps for pI-Code multiplexing:

1. Separation. Isolation of non-overlapping pI windows from each sample; Sample 1: peptides with pI < 5.0, Sample 2: peptides pI between 5.2-7.5, and Sample 3: peptides with pI > 7.8.
2. Pooling of all isolated pI ranges into one sample, followed by its LC MS/MS analysis.
3. Signal deconvolution: theoretical pI calculation of every peptide & assignment according to: pI < 5 to Sample 1, pI between 5.2 and 7.5 to Sample 2, pI > 7.8 to Sample 3.

Results

In isobaric labelling, the sample number information is encoded into a chemical tag, which is later decoded by MS [1]. Similarly, the isoelectric point value can be used as an information tag, since its value can be easily calculated, and because each protein produces tryptic peptides having pI values from acidic to neutral to basic, with empty valleys between them. In fact, it has been shown full human proteome coverage using peptides with pI ranges between 2.5 - 4.7 [2].

pI-Coding was developed as followed:

In-silico modelling. Several isoelectric focusing ranges were studied for optimum multiplexing. The results showed that two to three pI-channels provide at least 3 peptides for each protein.

pI-Coding. An isoelectric focusing device was developed for optimum collection of 3 pI areas in a way compatible with automation.

Validation. The results showed that pi-coding could be used for x2 or x3 multiplexing for the analysis of thousands of proteins in a highly parallelized manner for iTRAQ/TMT-based or label-free quantification.

Conclusions

The pi-Coding works with LFQ and TMT-reagents to double or triple their throughput. In LFQ, pi-coding is especially suitable for combining 3 biological replicates into a single sample (using 3 pi-Code channels). In the case of TMT-labeled samples, pi-Coding allows up to 20-30 sample multiplexing. In-silico modelling suggest that pi-Coding could be applied to different organisms. To our knowledge this is the highest level of multiplexing described in the literature.

Novel Aspect

pi-Coding works with LFQ and TMT-labeling to double or triple their throughput. In TMT-labeled samples, it allows up to 20-30 sample multiplexing.

References

1. Rauniyar N., and Yates J. *Journal of Proteome Research*, 13, 5293–5309 (2014).
2. Branca RM, et al. *Nature Methods*, 11, 59-62 (2014).

MP-157 / DIFFERENTIAL PROTEOMIC ANALYSIS OF ACTINIC KERATOSIS, BOWEN'S DISEASE AND CUTANEOUS SQUAMOUS CELL CARCINOMA BY SEQUENTIAL WINDOW ACQUISITION OF ALL THEORETICAL MASS SPECTRA (SWATH-MS)

Ali Azimi⁽¹⁾ - *Marina Ali*⁽¹⁾ - *Kimberley Kaufman*⁽²⁾ - *Pengyi Yang*⁽³⁾ - *Graham Mann*⁽⁴⁾ - *Pablo Fernandez-Penas*⁽¹⁾

Centre for Translational Skin Research, The University of Sydney, Westmead, Australia (1) - School of Life and Environmental Science, The University of Sydney, Camperdown, Australia (2) - School of Mathematics and Statistics, The University of Sydney, Camperdown, Australia (3) - Centre for Cancer Research, The University of Sydney, Westmead, Australia (4)

Keywords: Differential Proteomic Analysis of Actinic Keratosis, Bowen's Disease and Cutaneous Squamous Cell Carcinoma by Sequential Window Acquisition of All Theoretical Mass Spectra (SWATH-MS)

Introduction

Actinic keratosis (AK), Bowen's disease (BD) and cutaneous squamous cell carcinoma (cSCC) are heterogeneous skin tumours with high prevalence rates globally. This study reports comprehensive protein biomarkers with discriminatory capabilities towards these lesions.

Methods

Proteolytically digested peptides from formalin-fixed and paraffin embedded (FFPE) samples of normal skin, AK, BD, well differentiated cSCC (n=14), moderately differentiated cSCC (n=15) and poorly differentiated cSCC (n=15) were analysed in TripleTOF® 6600 using sequential window acquisition of all theoretical mass spectra (SWATH-MS) workflow. Peptide peaks were extracted using Skyline followed by their differential abundance analysis and bioinformatics analysis.

Results

We reliably identified and quantified 3574 proteins across 93 samples. Differential abundance analysis revealed that 15, 26 and 403 proteins were exclusively changed in cSCC, BD and AK samples compared to the normal skin, respectively. Also, 16 and 225 proteins were differentially abundant exclusively in BD and cSCC compared to AK samples, respectively. Ten proteins were differentially abundant in cSCC compared to BD samples. In addition, a total of 118 proteins were commonly changed in moderately-differentiated and poorly-differentiated cSCC samples compared to well-differentiated cSCC subtype while 55 and 809

proteins were differentially abundant exclusively in moderately and poorly differentiated cSCCs compared to well-differentiated cSCC samples, respectively. Majority of these proteins were involved in pathways associated with cell proliferation, apoptosis and DNA damage response which are crucial in carcinogenesis.

Conclusions

To our knowledge, this is the most comprehensive and complete SWATH-MS proteomic study in cSCC and its premalignant lesions. The identification of differentially abundant proteins and potentially important molecular pathways in this study form a rich resource for biomarker discovery in the future.

Novel Aspect

For the first time we successfully applied SWATH proteomics on FFPE skin samples. Our methodology opens invaluable opportunities to explore proteome of other skin conditions.

MP-158 / TAPE STRIPPED STRATUM CORNEUM SAMPLES PROVE TO BE SUITABLE FOR COMPREHENSIVE PROTEOMIC INVESTIGATION OF ACTINIC KERATOSIS

Ali Azimi⁽¹⁾ - *Marina Ali*⁽¹⁾ - *Kimberley Kaufman*⁽²⁾ - *Graham Mann*⁽³⁾ - *Pablo Fernandez-Penas*⁽¹⁾

Centre for Translational Skin Research, The University of Sydney, Westmead, Australia (1) - School of Life and Environmental Science, The University of Sydney, Camperdown (2) - Centre for Cancer Research, The University of Sydney, Westmead (3)

Keywords: tape stripping, mass-spectrometry based proteomics, actinic keratosis, noninvasive skin biopsy

Introduction

Actinic keratoses (AK) are frequent premalignant tumors that can be clinically difficult to differentiate from skin cancer. To definitely confirm or rule out the diagnosis, a biopsy is needed, but this is an invasive procedure that leads to increased morbidity. An easy, quick and reliable non-invasive alternative to biopsy is needed to definitively confirm the clinical diagnoses. In this paper, we evaluate Tape Stripping (TS) of stratum corneum (SC) for biomarker analysis of AK.

Methods

Lesional and non-lesional human SC samples were obtained by application of stripping tapes on the skin of five AK patients. Following sample preparation, protein digests were analyzed by liquid chromatography tandem mass-spectrometry (LC-MS/MS). Bioinformatics analysis were performed using Funrich, Ingenuity Pathway Analysis (IPA) and OncoPrint bioinformatics and analytical tools.

Results

Of the total 613 unique proteins identified, 477 overlapped with proteins identified in our proteomic analysis of formalin-fixed and paraffin-embedded (FFPE) AK samples. Additionally, 32 proteins were significantly increased, and four proteins decreased in AK samples compared to the normal skin ($p < 0.05$). In line with proteomic analysis of FFPE samples, IPA and Funrich analysis showed that differentially abundant proteins in the TS AK samples are implicated in PI3K/AKT and EGF signaling pathways. These findings were confirmed at the transcript level.

Conclusions

Tape stripped AK sample is suitable for biomarker analysis using mass spectrometry.

Novel Aspect

The application of this technique further could revolutionize management of keratinocytic skin tumors by reducing the need for traditional invasive biopsy.

MP-159 / A HIGH-RESOLUTION ORGANELLAR PROTEOMICS APPROACH TO STUDY PROTEIN SUBCELLULAR DISTRIBUTION AND DYNAMICS IN CELLS

Clara Steinbrueck - Johannes Jordan - Maciej Kocylowski - Astrid Kollwe - Wolfgang Bildl - Alexander Haupt - Bernd Fakler - Uwe Schulte

Institute of Physiology, Faculty of Medicine, University of Freiburg, Freiburg

Keywords: organellar proteomics, subcellular fractionation, relative quantification, hierarchical clustering, protein complexes

Introduction

In previous studies we have presented comprehensive interactome analysis using multi-epitope AP-MS ([1], [2]) and stoichiometric complexome analysis by csBN-MS ([3]) for membrane proteins. Notwithstanding, these techniques do not provide information on biogenesis of assemblies, subcellular distribution or spatiotemporal dynamics. We, therefore, took on extending/refining organellar proteomics workflows ([4], [5], [6], [7]).

Methods

Primary epithelial cells (RPTEC) were used as a model for endosomal recycling. After isotonic homogenization, hypotonic lysis and disruptive sonication multidimensional protein fractionation was applied. Fractions were analyzed by nano-LC-MS/MS. After label-free quantification, Euclidian-distance based hierarchical clustering and dimensionality reduction by t-distributed stochastic neighbour embedding (t-SNE) data were visualized in two dimensions.

Results

Key elements of the approach are subcellular fractionation (which is inherently limited), comprehensive and quantitative MS profiling, and unsupervised cluster analysis. High reproducibility and accuracy are required for optimal results. The higher the number of non-redundant fractions, the higher the precision of quantification. Nano-LC-MS/MS analysis of 80 trypsin-digested RPTEC samples (40 sub-fractions, each divided into an upper and a lower molecular weight range sample by SDS-PAGE) with an LTQ Orbitrap XL instrument yielded more than 3700 quantitative protein distribution profiles. Hierarchical clustering revealed more than 20 distinct entities representing classical cellular organelles or organelle substructures, populations of trafficking vesicles, multiprotein machineries and cytoskeletal networks as shown by successful mapping of 1074 prequalified marker proteins and validated functional annotations from UniProtKB/Swiss-Prot.

Conclusions

Organellar proteomics is a powerful method for comprehensive, robust and quantitative assessment of the subcellular localization of proteins with a broad range of applications: Identification and verification of subcellular markers, classification of interactomes, unbiased definition of subcellular proteomes, analysis of spatial dynamics of proteins and localization determinants, etc.

Novel Aspect

With multidimensional fractionation and advanced quantification, cluster analysis and visualization tools protein clusters were resolved with unprecedented resolution.

References

1. Müller C. S., Haupt A., Bildl W., Schindler J., Knaus H. G., Meissner M., Rammner B., Striessnig J., Flockerzi V., Fakler B., Schulte U., Proc. Natl. Acad. Sci. USA, 107, 14950-14957 (2010).
2. Schwenk J., Harmel N., Brechet A., Zolles G., Berkefeld H., Müller C. S., Bildl W., Baehrens D., Hüber B., Kulik A., Klöcker N., Schulte U., Fakler B., Neuron, 74, 621-633 (2012).

3. Müller C. S., Bildl W., Haupt A., Ellenrieder L., Becker T., Hunte C., Fakler B., Schulte U., Mol. Cell. Proteomics, 15, 669-681 (2016).
4. Foster L. J., de Hoog C. L., Zhang Y., Zhang Y., Xie X., Mootha V. K., Mann M., Cell, 125, 187-199 (2006).
5. Gilchrist A., Au C. E., Hiding J., Bell A. W., Fernandez-Rodriguez J., Lesimple S., Nagaya H., Roy L., Gosline S. J., Hallett M., Paiement J., Kearney R. E., Nilsson T., Bergeron J. J., Cell, 127, 1265-1281 (2006).
6. Itzhak D. N., Tyanova S., Cox J., Borner G. H., eLife, e16950 (2016).
7. Itzhak D. N., Davies C., Tyanova S., Mishra A., Williamson J., Antrobus R., Cox J., Weekes M. P., Borner G. H. H., Cell Rep., 20, 2706-2718 (2017).

MP-160 / EVALUATION OF PROTEOMICS QUANTITATIVE METHODS ON A RECONSTRUCTED SKIN MODEL

Sylvie Bourassa - Clarisse Gotti - Marie-Pier Scott-Boyer - Judith Marcoux - Isabelle Kelly - Arnaud Droit

CHU de Quebec Proteomics Platform, Laval University, Quebec, Canada

Keywords: proteomics quantification methods comparison softwares

Introduction

No consensus have been made whether label or label free proteomics quantitative methods are better for the quantitative analysis of complex mixtures. We have analyzed human reconstructed skin samples of 4 donors at different time using 10-plex TMT labeling, label free quantitation, DIA (Data Independent Acquisition) and SRM for a comprehensive analysis of the phenomena of skin reconstruction and comparison of proteomics quantitative techniques and softwares.

Methods

Epidermal skin reconstruction of 4 donors at 10 times (0 to 28 days) were used for the comparison. Peptides were labeled with 10-plex TMT and injected on the orbitrap fusion. The label-free experiment was injected on the orbitrap fusion and analysed by MaxQuant and Proteome Discoverer 2.2. 60 selected proteins were targeted by SRM on the 6500 QTRAP and analyzed by Skyline. DIA experiments were performed on the orbitrap fusion and analysed by Scaffold DIA or Skyline. .

Results

4338 proteins were identified by using 5 different 10-plex TMT experiments (14 fractions of 2 hours) with a common reference pool. 1821 quantifiable proteins were common between the 5 experiments and 48 variant proteins were obtained. For the Label-Free experiment, individual samples were run for 5 hours identifying 3619 proteins, from which 2749 were quantifiable. 393 proteins were considered variant with MaxQuant and the results will be compared with PD 2.2 software. 55 of the selected proteins for the SRM experiment were quantifiable and 13 proteins were variants. DIA experiments (50 windows of 10 amu) using Scaffold DIA and Skyline will be compared with the 3 other methods. 11 proteins were commonly statistically variants between the 4 different techniques at specific time points and 8 of them behave exactly the same along the time course for the quantitative techniques tested.

Conclusions

For this project, TMT gave the highest number of protein identifications, Label-Free analysis gave the highest number of quantifiable proteins and DIA gave the best results for the reproducibility. Label free analysis using precursor intensity was the most valuable method despite the longest time of analysis. We found that the use of Scaffold DIA is promising for the analysis of DIA data.

Novel Aspect

Different quantitative proteomics methods and the use of a new released software for DIA allows better understanding and confidence in the mechanism of skin reconstruction.

MP-161 / NEW PHARMACOLOGICAL TARGETS FOR CYSTIC FIBROSIS TREATMENT FROM EXTENSIVE PROTEOMIC PROFILING OF Δ F508-CFTR EXPRESSING CELLS

Clarissa Braccia⁽¹⁾ - *Valeria Tomati*⁽²⁾ - *Nicoletta Pedemonte*⁽²⁾ - *Andrea Armirotti*⁽³⁾

Fondazione Istituto Italiano di Tecnologia, D3 Pharmachemistry, Genova (1) - Istituto Giannina Gaslini, U. O. C. Genetica Medica, Genova (2) - Fondazione Istituto Italiano di Tecnologia, Analytical Chemistry Facility, Genova (3)

Keywords: Cystic Fibrosis; bronchial epithelial cells; SWATH ion library

Introduction

Cystic fibrosis (CF) is a genetic disorder, caused by mutations in the CFTR gene [1]. Deletion of phenylalanine 508 (F508del) is the most frequent mutation. Dr. Pedemonte's group has selected molecules with a role in F508del-CFTR processing and degradation. After their silencing, an increase of CFTR at the cell membrane was measured [2], [3]. Using SWATH workflow [4], it is possible to analyze changes at proteome level. The first step is the creation of an ion library.

Methods

We characterized the proteome of CFBE410-, the in-vitro model for cystic fibrosis pharmacological research [5]. CFBE410- cells were cultured and transfected as already described [2]. From cell lysates, we performed both an in-gel digestion and an in-solution digestion. The in-solution digestion was performed on both Proteominer [6] wash and eluted fractions and the resulting peptides were then fractionated offline (8 fractions) with a high pH/low pH strategy [7].

Results

Raw data were analyzed with ProteinPilot software using the Paragon algorithm [8]. Spectra were searched against the Homo sapiens reference proteome, reporting 71599 proteins. Search was performed against both target and decoy databases. Carbamidomethylation of cysteine (CAM) was set as fixed modification and only four variable modifications were allowed, maximum one per peptide (methionine oxidation, asparagine and glutamine deamidation and conversion of glutamine and glutamate to pyroglutamic acid). This search, done using the FDR calculation protocol described by Tang in 2008 [9], returned 5935 proteins and 72031 distinct peptides at a FDR of 1%.

To ensure high quality to our library, we imported in the library builder only MRM data of proteins with a maximum global FDR of 1%. Moreover, we removed the peptides shared by more than one protein. After these filtering steps, we produced an ion library with the spectral information on 5627 proteins.

Conclusions

Thanks to this ion library, we were able to quantify more than 2000 proteins in a total lysate from a single 2 hours LC-MS/MS run.

The ion library we generated will enable us to detect changes at protein level in response to the silencing of molecules known to be involved in CFTR processing.

Novel Aspect

We built the first SWATH ion library dedicated to the most used cell model for CF research. Unlike the pan-human ion library [10], our library contains CFTR and some of CFTR interactors.

References

1. Kumar, S., A. Tana, and A. Shankar, Cystic fibrosis-what are the prospects for a cure? *Eur J Intern Med*, 2014. 25(9): p. 803-7.
2. Tomati, V., et al., Genetic Inhibition Of The Ubiquitin Ligase Rnf5 Attenuates Phenotypes Associated To F508del Cystic Fibrosis Mutation. *Sci Rep*, 2015. 5: p. 12138.

3. Tomati, V., et al., High-throughput screening identifies FAU protein as a regulator of mutant cystic fibrosis transmembrane conductance regulator channel. *J Biol Chem*, 2017.
4. Anjo, S. I., C. Santa, and B. Manadas, SWATH-MS as a tool for biomarker discovery: From basic research to clinical applications. *Proteomics*, 2017. 17(3-4).
5. Avramescu, R. G., et al., Mutation-specific downregulation of CFTR2 variants by gating potentiators. *Human Molecular Genetics*, 2017. 26(24): p. 4873-4885.
6. Hartwig, S., et al., Combinatorial hexapeptide ligand libraries (ProteoMiner™): An innovative fractionation tool for differential quantitative clinical proteomics. *Archives of Physiology and Biochemistry*, 2009. 115(3): p. 155-160.
7. Wang, H., et al., An off-line high pH reversed-phase fractionation and nano-liquid chromatography-mass spectrometry method for global proteomic profiling of cell lines. *J Chromatogr B Analyt Technol Biomed Life Sci*, 2015. 974: p. 90-5.
8. Shilov, I. V., et al., The Paragon Algorithm, a next generation search engine that uses sequence temperature values and feature probabilities to identify peptides from tandem mass spectra. *Mol Cell Proteomics*, 2007. 6(9): p. 1638-55.
9. Tang, W. H., I. V. Shilov, and S. L. Seymour, Nonlinear fitting method for determining local false discovery rates from decoy database searches. *J Proteome Res*, 2008. 7(9): p. 3661-7.
10. Rosenberger, G., et al., A repository of assays to quantify 10,000 human proteins by SWATH-MS. *Sci Data*, 2014.1: p. 140031.

MP-162 / STRAIGHT FORWARD PHOSPHOPROTEOMIC ANALYSIS PIPELINE IN CHO CELLS AS MAMMALIAN NON-MODEL ORGANISM

Louise Brachtvogel⁽¹⁾ - *Antonia Pries*⁽¹⁾ - *Stefan Walter*⁽²⁾ - *Thomas Noll*⁽¹⁾ - *Raimund Hoffrogge*⁽¹⁾

Cell Culture Technology, Bielefeld University, Bielefeld (1) - Faculty of Biology, University of Osnabrück, Osnabrück (2)

Keywords: CHO, MaxQuant, LFQ, SILAC, Signaling

Introduction

Quantitative (phospho)proteomics is often used to elucidate cellular responses on media supplements [1]. Especially for non-model species, certain aspects besides LC-MS set up, like sample preparation, databases and evaluation workflow strongly impact the outcome of signaling studies. Adapted workflows then lead to impacted phosphorylation sites and signaling pathways. Here we outline critical steps based on our experience with chinese hamster ovary (CHO) cell analysis.

Methods

CHO cells were cultivated in chemically defined media. For SILAC experiments, media are supplemented with either lLys/lArg or hLys/hArg. After protein digestion, phosphopeptides were enriched with different strategies. Measurements were performed by nLC-ESI Orbitrap MS and data was evaluated utilizing two different pipelines based on MaxQuant/Perseus and Proteome Discoverer (PD). The obtained SILAC data was also quantified by LFQ for cell line comparison.

Results

Sample preparation workflow

An optimized lysis buffer solubilizes membrane proteins and allows a direct digest. The resulting new 6-step workflow decreases the variation of identified phosphorylation sites between techn. replicates.

Database searches

Combination of two databases (well annotated mouse; incomplete *C. griseus* for CHO) leads to the highest identification rate compared to individual DB. However, the mix of identifiers for both species complicates

software evaluation. Merging of data fore. g. pathway mapping hence is to great extend manually. An intermediate step was necessary for mapping different protein identifiers to homologous ones.

Quantification Methods

The SILAC-based data analysis yields in substantially more identification of modified peptides then LFQ-based approaches, where only 75 % were identified. However, LFQ offers to describe cell-specific ground phosphorylation status and protein expression comparisons. As the current version of PD supports LFQ, we summarize pros and cons in comparison to MaxQuant.

Conclusions

Several aspects besides the LC-MS set up showed to have strong impact on signaling analysis of non-model organisms such as CHO. We highlight the importance and limitations of sample preparation steps, referred protein databases and quantification methods and propose applicable solutions.

With this improvements, based on complementary LFQ + SILAC phosphoproteomic data, a comprehensive view on CHO cell signaling in dependence on media supplements was possible.

Novel Aspect

Optimized procedures allowed us to reach a high coverage of less-abundant phosphoproteins in CHO cells, which results in a full view on signaling events in this non-model organism.

References

1. Müller B., Heinrich C., Jabs W., Kaspar-Schönefeld S., Schmidt A., Rodrigues de Carvalho N., Albaum S. P., Baessmann C., Noll T., Hoffrogge R., J Biotechnol., 257:87-98 (2017).

MP-163 / COMPARATIVE PROTEOMICS STUDY FOR ELUCIDATION OF POTENTIAL MECHANISMS UNDERLYING YEAST RESISTANCE TO ALCOHOLS

Julia Bubis⁽¹⁾ - *Vladimir Gorshkov*⁽²⁾ - *Frank Kjeldsen*⁽²⁾ - *Mikhail Gorshkov*⁽¹⁾ - *Irina Tarasova*⁽¹⁾ - *Daria Spasskaya*⁽³⁾ - *Dmitry Karpov*⁽³⁾

Talrose Institute for Energy Problems of Chemical Physics, Moscow (1) - Department of Biochemistry and Molecular Biology, University of Southern Denmark, Odense (2) - Engelhardt Institute of Molecular Biology, Moscow (3)

Keywords: shotgun proteomics, comparative proteomics, *S. cerevisiae*, alcohol stress

Introduction

Distilled spirits production using *S. cerevisiae* requires understanding the mechanisms of yeast cell resistance to alcohols. Mutations in genes of ubiquitin-proteasome proteolytic system, such as RPN4, may result in appearance of the strains exhibiting hyperresistance [1]. Indeed, RPN4 controls expression of proteolytic genes involved in metabolic pathways [2]. In this work we elucidate the RPN4-dependent mechanisms of yeast cell resistance to alcohol.

Methods

Three *S. cerevisiae* strains: wild-type (WT) BY4742, RPN4-Δ (with RPN4 deletion), and YPL mutant with decreased proteasome activity, were grown up to log phase in stressed (7% ethanol) and normal conditions. Samples were subjected to LC-MS/MS analysis followed by protein identification, quantification (X!Tandem/MPscore), statistical evaluation of changes in protein expression, and pathways analysis (STRING).

Results

PCA results have shown the protein grouping according to strains, biological and technical replicates. No

technical bias was found for analyzed data. Next, protein levels were compared in pairs: no stress vs. ethanol stressed for each strain (WT, YPL, RPN4-Δ). Proteins that passed the corrected (Benjamini–Hochberg FDR) t-test's p-value of 0.05 were defined as differentially expressed. We found the following pathways affected significantly by the ethanol effect on the cells: protein folding and unfolding, rRNA modification, as well as carbohydrate and nucleotide metabolism. Thus, the ethanol elicits complex cellular response that involves known players like chaperones and metabolic enzymes, as well as the new ones like enzymes of nucleotide metabolism and rRNA modification. Some of these genes are RPN4-dependent. Detailed analysis of differentially expressed proteins allowed identifying the candidates that determine the phenotype of the cells' alcohol resistance.

Conclusions

Comparative proteome analysis of three different strains of *S. cerevisiae* revealed differences in yeast response to ethanol at proteome level and suggested gene candidates that are potentially responsible for developing alcohol resistance in the cells. These genes can be targets for further studies using methods of molecular genetics.

Novel Aspect

The importance of RPN4-regulated non-proteasomal genes as novel factors involved in cellular response to ethanol was revealed for the first time using LC-MS/MS.

The work was supported by the Russian Foundation for Basic Research (grant #18-04-00692) and Russian Science Foundation (project #14-14-00971).

References

1. Gonzalez-Ramos D., et al. *Biotechnology for Biofuels*, 6:48 (2013).
2. Mannhaupt, G. et al. *FEBS Lett.* 450: 27-34 (1999).

MP-164 / QUANTITATIVE LABEL-FREE PROTEOMICS IN A NOVEL CELLULAR KNOCK OUT MODEL OF METHYLMALONIC ACIDEMIA

Michele Costanzo⁽¹⁾ - ***Marianna Caterino***⁽¹⁾ - ***Armando Cevenini***⁽¹⁾ - ***Ida Chiara Guerrero***⁽²⁾ - ***Margherita Ruoppolo***⁽¹⁾

University of Naples "Federico II", Department of Molecular Medicine and Medical Biotechnology, Naples (1) - Plateforme Protéomique 3P5-Necker, INSERM U24, Structure Fédérative de Recherche Necker and Université Paris Descartes, Paris (2)

Keywords: methylmalonicacidemia, quantitative proteomics, mass spectrometry, bioinformatics, CRISPR-CAS9

Introduction

Methylmalonic Acidemias (MMAs) are severe autosomal recessive inborn errors of metabolism caused by deficiency of the mitochondrial methylmalonyl-CoA mutase (MUT) enzyme, which converts methylmalonyl-CoA into succinyl-CoA, a Krebs cycle intermediate. As consequence, methylmalonic acid accumulates in body fluids as biomarker of the disease. Long-term complications can include cognitive deficit, chronic kidney disease, liver and kidney failure.

Methods

A stable MUT-knock out (KO) cell model was generated by using CRISPR-CAS9 genome editing technology on the embryonic kidney (HEK) 293 cell line. Methylmalonic acid was measured in KO cells by liquid chromatography - tandem mass spectrometry (LC-MS/MS) in the multiple reaction monitoring (MRM)

mode. The proteomes of 4 biological MUT-KO and 4 wild type (WT) cell samples were compared using a quantitative label-free approach by LC-MS/MS, according to MaxQuant software.

Results

The goodness of this new cellular model, in which MUT protein was completely knocked out, was tested measuring the levels of methylmalonic acid inside the cells. By contrast, it has not been detected in WT cells, as expected.

Bioinformatics analysis was carried out downstream the proteomic experiment in order to select the significant differentially expressed proteins between the two conditions, among those identified in the 8 samples. The analysis showed a number of 178 differentially expressed proteins, which passed all the statistic tests (results from Perseus platform). Applying a more stringent parameter (the identification of each protein with at least 5 peptides in at least 3 of 4 samples per condition), the number of significant proteins was further reduced to 69, of which 39 down-regulated and 30 up-regulated in the MUT-KO condition.

Conclusions

The selected proteins are globally involved in energy metabolism and oxidation-reduction processes and most of them have mitochondrial localization, supporting the idea of mitochondrial impairment in MMAs. The rescue of MUT protein will be performed in KO cells in order to specify if the deregulated proteins are strictly linked to MUT absence and highlight unknown molecular mechanisms underlying the cellular damage and tissue alterations in MMAs.

Novel Aspect

This is the first cellular model of methylmalonicacidemia generated by using CRISPR-CAS9 genome editing technology.

MP-165 / INACTIVATED POLIOVIRUS VACCINE ANTIGEN DETERMINATION: A PROTEOMICS APPROACH

Marybeth Creskey - Diane Viel - Tong Wu⁽³⁾ - Terry Cyr

Health Canada, Biologics and Genetic Therapies Directorate, Ottawa

Keywords: vaccine, quantitation, proteomics, bottom-up, antigen

Introduction

Polio (poliomyelitis) once a widespread, occasionally paralytic diseases almost eradicated. Key to this success has been broad vaccination programs with oral poliovirus vaccine (OPV) but the endgame requires a switch to inactivated vaccines (IPV)(1). IPV potency is assayed in vivo(immunogenicity)(2) or in vitro (antibody based)(3) by relative comparison to international references; improved IPV testing methods are desired(4).

Methods

Reference Antigens (RA), Inactivated Polio Vaccine (IPV), monovalent bulks, and combination vaccines (DTaP-IPV) were disulfide cleaved and capped then digested with trypsin. The peptides were separated on a C18 column connected to Orbitrap Fusion MS collecting MS/MS in data dependent mode with HCD fragmentation. The data was analyzed with Proteome Discoverer 2.1 by searching a custom database and quantification of proteins based on top 3 precursor intensity(5).

Results

We compared numerous polio standards and vaccines using proteomics and label free quantitation. The database search required generation of a custom database consisting of individual polio capsid proteins

and proteins derived from the other components of the combination vaccine. A recent RA with established DU/ml values was used to prepare a calibration curve and the equation relating precursor intensities to DU was determined. The calculated Dus in the samples were then compared to the results obtained by the manufacturers by antibody based methods. The DU values obtained by this method are higher than labelled amounts for the combination vaccines, IPV, and reference antigens while the bulks were lower than the labelled content. An independent preparation of the calibration curve RA has calculated DU values close to expected. Identification and quantification of polio, diphtheria, pertussis, tetanus proteins in combination vaccines is possible.

Conclusions

The label free quantitation of Reference Antigens using optimized MS methods and bioinformatics software provides results that are closely aligned with those from antibody based methods. The advantage of the proteomics approach is that unexpected or adventitious proteins in the sample are also identified and quantified without prior knowledge of their presence. The disadvantage of the proteomics approach is the lack of proof regarding the protein folding.

Novel Aspect

A rapid proteomics and bioinformatics approach has been developed for polio vaccine antigen quantitation.

References

1. Zaffran M, McGovern M, Hossaini R, Martin R, Wenger J. The polio endgame: securing a world free of all polioviruses. *Lancet*. 2018;391(10115):11-3.
2. Van Steenis G, Van Wezel AL, Sekhuis VM. Potency testing of killed polio vaccine in rats. *Dev Biol Stand*. 1980;47:119-28.
3. Souvras M, Montagnon B, Fanget B, van Wezel AL, Hazendonk AG. Direct enzyme linked immunosorbent assay (ELISA) for quantification of poliomyelitis virus D-antigen. *Dev Biol Stand*. 1980;46:197-202.
4. Baca-Estrada M, Griffiths E. Regulation and standardization of IPV and IPV combination vaccines. *Biologicals*. 2006;34:159-61.
5. Creskey MC, Li C, Wang J, Girard M, Lorbetskie B, Gravel C, et al. Simultaneous quantification of the viral antigens hemagglutinin and neuraminidase in influenza vaccines by LC-MS E. *Vaccine*. 2012;30:4762-70.

MP-166 / TMT-BASED QUANTITATIVE APPROACH DISCLOSES MONACOLIN K-INDUCED PROTEOME ALTERATIONS IN TRIPLE NEGATIVE BREAST CANCER CELLS

Federica Del Gaudio⁽¹⁾ - **Maria Chiara Monti**⁽¹⁾ - **S. Rontogianni**⁽¹⁾ - **Kelly Steacker**⁽²⁾ - **Ida Chiara Guerrero**⁽³⁾ - **Raffaele Riccio**⁽¹⁾ - **Albert J. R. Heck**⁽⁴⁾ - **Maarten A. F. Altelaar**⁽⁴⁾

Department of Pharmacy, University of Salerno, Via Giovanni Paolo II, 132, 84084 Fisciano (SA), Italy (1) - Biomolecular Mass Spectrometry and Proteomics, Bijvoet Center for Biomolecular Research and Utrecht Institute for Pharmaceutical Sciences, Utrecht University, Padualaan 8, 3584 CH Utrecht, The Netherlands. (2) - Paris Descartes-Sorbonne Paris Cité University, Paris, France (3) - Netherlands Proteomics Centre, Padualaan 8, 3584 CH Utrecht, The Netherlands (4)

Keywords: Quantitative Mass Spectrometry, TMT-labeling, Natural compounds, Nutraceuticals, Triple Negative Breast cancer (TNBC).

Introduction

MS-based Proteomics on Natural Compounds.

MS-based proteomics represents a cleaver tool to disclose unclear bioactivities of small molecules such as natural compounds in nutraceuticals [1]. The latter are rising as significant building block in prevention or/and treatment of numerous diseases [2].

Monacolin K (MNK), abundant in red yeast rice extract, is used to lower hematic concentration of cholesterol. Recent scientific evidences link statins to antiproliferative and apoptotic effects in a wide panel of cancers [3].

Methods

TMT-based quantitative proteomics.

A kinetic proliferation assay was performed using IncuCyte ZOOM® system to find out the IC50 of MNK to set up proteomics experiments.

Proteomes obtained upon different treatments of MDA-MB231 cells with MNK were digested and pooled together, after former TMT 10-plex™ labeling of the related peptide mixtures [4][5].

The obtained multiplexed pool was prepared for the UPLC-MS/MS analysis of proteomes and phosphoproteomes. Bioinformatic analysis were carried out by MaxQuant and Perseus [6][7].

Results

Monacolin K affects proteome and phosphoproteome of MDA-MB231.

MDA-MB231 cells were treated with different concentrations of MNK and growth kinetic curves unraveled an IC50= 4.02±1.23µM.

This time-course experiment was the groundwork to drug treatments for the proteomics approach: MDA-MB231 cells were treated with MNK at 5 µM or 0.5 µM and harvested after 30', 6h or 24h. Control experiments without MNK were harvested in parallel. Then, cells were lysate and proteomes were digested. The obtained tryptic mixtures were chemically labeled with amino-reactive tandem mass isobaric tags (TMT 10-plex™). Labeled peptides were pooled together and the multiplexed sample was simplified upon HpH-fractionation. An aliquot of each fraction was used to obtain the phosphoproteome upon Fe³⁺-NTA affinity chromatography. NanoUPLC-MS2 experiments were carried out both for proteomes and phosphoproteomes. Raw MS data were analyzed with MaxQuant software to identify and quantify proteins. Statistical analysis of MQ results was performed by Perseus. When MDA-MB231 are treated with MNK 5µM, some alterations in their proteome and phosphoproteome have been found.

Conclusions

Monacolin K regulates proliferative pathways in TNBC cells.

MNK has an antiproliferative effect on MDA-MB231, as shown by growth kinetic curves. MS-based proteomics unrevealed the potentiality of this bioactive metabolite, abundant in red yeast rice, in modulating expression and phosphorylation of proteins involved in cellular pathways responsible for cellular proliferation. These findings configure Monacolin K as an interesting metabolite in adjuvant therapies for the treatments of triple negative breast cancer.

Novel Aspect

Exploiting TMT-10plex™, a quantitative MS-based proteomics approach was used to investigate antiproliferative effect of the natural compound Monacolin K on MDA-MD231 cells.

References

1. Rix, U., Superti-furga, G., Nat. Chem. Biol., 5, 616-624 (2009).
2. Aronson, J. K., Br. J. Clin. Pharmacol., 83, 8-19 (2017).
3. Ahern, T. P., Damkier, P., Cronin-Fenton, D. P., Thomas, C. T., Ahern, P., Lancet Oncol, 151016, 461-468 (2014).
4. Thompson, A., Schalfer, J., Kuhn, K., Kienle, S., Schwarz, J., Schmidt, G., Neumann, T., Hamon C., Anal. Chem., 75, 1895-1904 (2003).
5. Altelaar, A. F. M. et al., J. Proteomics, 88, 14-26 (2013).
6. Cox, J., Mann, M., Nat. Biotechnol., 26, 1367-72 (2008).
7. Tyanova, S., Temu, T., Sinitcyn, P., Carlson, A., Hein, M., Geiger, T., Mann, M., Cox, J., Nature Methods, 13,731-740 (2016).

MP-167 / HIGH THROUGHPUT SIGNALING PATHWAY MONITORING BASED ON FAST LC-PRM ANALYSES

Sebastien Gallien^(1,2) - Aaron Gajadhar⁽³⁾ - Bhavin Patel⁽³⁾ - Tabiwang N. Arrey⁽⁴⁾ - Dave Sarracino⁽²⁾ - Yue Xuan^(2,4) - Emily Chen⁽²⁾

Thermo Fisher Scientific, Paris, France; (1) - Thermo Fisher Scientific, Precision Medicine Science Center, Cambridge, MA, USA; (2) - Thermo Fisher Scientific, San Jose, CA, USA (3) - Thermo Fisher Scientific, Bremen, Germany (4)

Keywords: Parallel reaction monitoring; signaling pathway; high throughput

Introduction

Targeted analyses based on HRAM parallel reaction monitoring (PRM) measurements have opened new opportunities in quantitative proteomics. Here PRM was combined with fast capillary-flow LC separation to accelerate the execution of analyses. The setups developed turned out to be ideally suited to high throughput monitoring of signaling pathways in the context of clinical research applications, typically encompassing thousands of samples.

Methods

The analyses were performed on a Thermo Scientific™ Q Exactive™ HF mass spectrometer operated with several PRM-based acquisition schemes (using instrument programming interface in some cases). Chromatographic separation were carried out using a Thermo Scientific™ UltiMate™ 3000 RSLC system equipped for capillary flow and a Evosep-One system (Evosep). Various gradient lengths and MS parameter settings were employed to analyze samples of various complexities.

Results

Fast LC methods have been designed with the objective to increase the analytical throughput, therefore requiring both higher flow rate (typically 2-3 $\mu\text{L}/\text{min}$) and gradient length reduction, while maintaining acceptable duty cycle. Different experiments have been developed for conventional PRM, dynamic retention time-PRM (dRT-PRM), and internal standard triggered-PRM (IS-PRM) analyses of ten to a few hundred endogenous peptides. Another IS-PRM method variant, relying on the concomitant measurement of each pair of internal standard and endogenous peptides in a single scan, was explored to further shorten analysis time. The final methods differentiated themselves by their unique combination of acquisition scheme, experiment scale, gradient length, and MS parameter settings. They were applied to the analyses of samples of high complexity, *e.g.*, digests of cell lines, optionally mixed with Yeast digest, and samples of moderate complexity obtained through multiplex immunoprecipitation targeting proteins of AKT/mTOR pathway.

Conclusions

A detailed comparative evaluation of the quantification performance (*e.g.*, quantification precision, accuracy, and sensitivity) obtained for various experimental configuration was carried out. This in-depth investigation enabled the establishment of a reference baseline demonstrating sensitive quantification of a few dozens of peptides in 100 samples within one day under high acquisition efficiency (typically through dRT-PRM and especially IS-PRM mode).

Novel Aspect

Sensitive PRM-based signaling pathway monitoring with the unique combination of analytical throughput and experiment scale required for clinical research applications.

MP-168 / QUANTITATION OF AMYLASE/TRYPsin-INHIBITORS IN VARIOUS WHEAT SPECIES

Sabrina Geisslitz⁽¹⁾ - Katharina Scherf⁽¹⁾ - Peter Koehler⁽²⁾

Leibniz-Institute for Food Systems Biology at the Technical University of Munich, Freising, Germany (1) - Biotask AG, Technical University of Munich, Esslingen am Neckar (2)

Keywords: α -amylase/trypsin-inhibitors, non-celiac gluten sensitivity, ancient and modern wheat species, stable isotope dilution assay, intensity based absolute quantitation

Introduction

The “ancient” wheat species einkorn, emmer and spelt are claimed to have sensory and health benefits compared to “modern” bread and durum wheats [1]. To date, it is not known if ancient wheat species contain lower contents of α -amylase/trypsin-inhibitors (ATI), which have been associated with intestinal inflammation and symptoms typical of non-celiac gluten sensitivity (NCGS) [2]. Therefore, the aim of this study was to develop a method to quantitate ATI in wheat.

Methods

ATI were extracted from eight wholegrain flours of the five wheat species with salt solution. Proteins were reduced, alkylated and digested with trypsin. Peptides were analyzed by LC-MS/MS and five ATI-types were quantitated by stable isotope dilution assays (SIDA). For the validation of trueness, the results of SIDA were compared with an independent LC-MS/MS system and data evaluation by intensity based absolute quantitation algorithm (iBAQ).

Results

Precision ($< 8\%$), recovery (80 – 120 %) and limits of detection ($< 5\text{ mg/kg}$) showed the good performance of the method. The accuracy of SIDA was confirmed by iBAQ, which yielded comparable results. In contrast to the other wheat species, only three of eight einkorn cultivars contained ATI in very low concentrations ($0.2 \pm 0.1\text{ mg/g}$). On average, all spelt ($4.7 \pm 0.5\text{ mg/g}$) and emmer cultivars ($5.0 \pm 0.5\text{ mg/g}$) had higher ATI contents than bread wheat ($3.7 \pm 0.3\text{ mg/g}$) and durum wheat ($4.2 \pm 0.9\text{ mg/g}$). The ATI content was not correlated with the total crude protein content ($r = 0.01$). Finally, differences in the distribution of ATI-types were found. The percentages of 0.19 and 0.28 ATI to CM-type ATI was balanced (52%/48 %) in bread wheat and spelt, whereas CM-types were predominant (25%/75 %) in durum wheat and emmer.

Conclusions

We developed a novel sensitive LC-MS/MS method based on SIDA for the quantitation of ATI in ancient and modern wheat species. Applying the new method, the assumption that ancient wheats contain less ATI than modern ones was disproven for spelt and emmer. The results indicated that different wheat species can be distinguished by the distribution of ATI. Comparison with iBAQ showed that one internal standard was sufficient for SIDA.

Novel Aspect

The ancient wheat species spelt and emmer contain more ATI than modern ones. Therefore it is not confirmed that ancient wheat species might be more suitable for NCGS patients.

References

1. Longin C. F. H., Ziegler J., Schweiggert R. et al., *Crop Science*, 10, 1-10 (2016)
2. Zevallos V. F., Raker V., Tenzer S. et al., *Gastroenterology*, 152, 1100-1113 (2017)

MP-169 / PEPTIDE MIGRATION PROFILES: A TOOL TO TRACE PROTEOLYTIC PROCESSES - FACTOR XII CONTACT SYSTEM ACTIVATION

Laura Heikau - Hartmut Schlüter

University Medical Center Hamburg-Eppendorf, Mass Spectrometric Proteomics, Hamburg

Keywords: Proteolysis, Factor XII, coagulation, LC-MS/MS, SDS-page

Introduction

Many signaling pathways are activated when proenzymes are proteolytically cleaved releasing the active enzyme (*e.g.* auto-activation of Factor XII by contact with negatively charged surfaces, starting pathways of coagulation and inflammation). In this work, we use migration profiles [1] for differential analysis of Factor XII contact activation by dextran sulfate (DXS) and poly-phosphate (polyP) and downstream effects.

Methods

Human Plasma was activated using DXS or polyP. An immunoprecipitation with anti-FXII antibody was performed on the activated plasma samples and a control sample. Eluates were loaded onto a 10% SDS-page. Each of the obtained lanes was cut into 19 bands for in gel tryptic digestion and LC-MS/MS analysis. Migration profiles were obtained by plotting the relative intensity of the identified peptides against the corresponding empirical molecular weight.

Results

The obtained migration profiles show a distinguished pattern for the FXII peptides in plasma activated with PolyP or DXS and in the control. When plasma is activated by DXS, it shows a clear separation of the heavy chain from the light chain, this separation is not observed in the case plasma was activated by polyP and in the control. In both cases of activated plasma, the downstream proteolytic processes are activated. Compared to the control, the activation of plasma kallikrein (PK) and release of bradykinin from high molecular weight kininogen (HK) is clearly demonstrated. These results corroborate previous findings by Engels et al. that suggest activation of FXII by polyP to be due to conformational changes, and not by cleavage at the activation site (R353-V354) [2].

Conclusions

Migration profiles of all identified proteins from differentially activated plasma allowed a comprehensive study of the activated proteolytic pathways. It's shown that the use of different activators has an impact on the mechanism of FXII activation. These findings might provide an explanation for FXII activation by polyP resulting in coagulation and activation of the kallikrein-kinin system, whereas by DXS only activates the kallikrein-kinin system [3].

Novel Aspect

The method presented is a useful tool to understand the differential activation of FXII and can be applied on a wide array of fields where proteolysis is an issue.

References

1. Kwiatkowski M. ; J Proteomics. 134, 193-202 (2016)
2. Engels R. . ; Journal of Thrombosis and Haemostasis, 12, 1513-1522 (2014)
3. Weidmann, H. ; BBA - Molecular Cell Research, 1864, 2118-2127 (2017)

MP-170 / MYELIN PROTEOMICS 2.0: APPLICATION TO MOUSE MODELS FOR MYELIN-RELATED DISEASES

Olaf Jahn⁽¹⁾ - **Siems Sophie**⁽²⁾ - **Dörte Hesse**⁽¹⁾ - **Stefan Tenzer**⁽³⁾ - **Hauke Werner**⁽²⁾

Max-Planck-Institute of Experimental Medicine, Proteomics Group, Göttingen (1) - Max-Planck-Institute of Experimental Medicine, Department of Neurogenetics, Göttingen (2) - University of Mainz, Institute of Immunology, Mainz (3)

Keywords: myelin, FASP, automation, quantification, neuropathies

Introduction

Myelin, a plasma membrane specialization of glial cells, facilitates rapid signaling in the nervous system. A key tool to approach failing glia-axonal support in neurological disorders is the proteome analysis of myelin, but membrane protein content and dynamic range pose technical challenges (reviewed in [1,2]). We show

advances in label-free quantification of myelin proteins and applications to a mouse model for a heritable peripheral neuropathy.

Methods

Myelin was solubilized with ASB-14 and subjected to FASP-based digestion automated on a liquid handling workstation (Tecan Freedom EVO 150). Data were acquired on a QTOF LC-MS-system (Waters Synapt G2-S) using three different data-independent acquisition (DIA) strategies: MSe, without ion mobility separation (IMS); UDMSe, with IMS; DRE-UDMSe, with IMS and dynamic range enhancement. Label-free protein quantification was performed with the IsoQuant software.

Results

As a proof-of-principle for robustness and accuracy of quantification, we first compared CNS myelin from wild type mice (CnpWT/WT) with myelin from mutant mice heterozygous for the CNP gene (CnpWT/null). In three independent experiments, each comprising 24 MSe runs (three biological replicates per condition, duplicate digestion and injection), we determined the ratio of the abundance of CNPCnpWT/WT:CnpWT/null as 1.98:1 (CV < 10%), exactly matching the ratio of 2:1 as expected from Cnp gene heterozygosity. In UDMSe runs, however, levels of high-abundant myelin proteins including CNP were misquantified due to the compression of the dynamic range caused by IMS of peptides. This issue was addressed by DRE-UDMSe, a novel data acquisition method in which a deflection lens is used to cycle between full and reduced ion transmission during one scan. All three data acquisition methods were applied to PNS myelin from mice deficient for periaxin (Prxnull/null), a model for a demyelinating neuropathy of the Charcot-Marie-Tooth type (CMT-4F).

Conclusions

While MSe and UDMSe are well suited for accurate quantification of highly abundant proteins and deep proteome coverage, respectively, DRE-UDMSe provides a compromise between dynamic range and identification rate needed for routine quantification of myelin proteins. Analysis of Prxnull/null mice systematically revealed altered abundance of potentially disease-relevant PNS myelin proteins and thus allows insights into myelin-related pathology in this model.

Novel Aspect

Development and automation of a customized FASP procedure; improved label-free quantification of myelin proteomes; routine proteomics phenotyping of mutant mice

References

1. Jahn O., Tenzer S., Werner H. B., Mol Neurobiol, 40, 55 (2009).
2. de Monasterio-Schrader P., Jahn O., Tenzer S., Wichert S. P., Patzig J., Werner H. B., Cell Mol Life Sci, 69, 2879 (2012).

MP-171 / SIMPLIFYING THE USE OF ION LIBRARIES DURING DATA PROCESSING OF DATA INDEPENDENT ACQUISITION PROTEOMICS DATA

Matthew Huebsch⁽¹⁾ - **Dietmar Waidelich**⁽²⁾ - **Adam Lau**⁽¹⁾ - **Kathleen Lewis**⁽³⁾ - **Nick Morrice**⁽⁴⁾ - **Jörg Dojahn**⁽²⁾ - **Dietrich Merkel**⁽²⁾ - **Rebekah Sayers**⁽⁴⁾ - **Zuzana Demianova**⁽²⁾ - **Sara Ahadi**⁽⁵⁾ - **Christie Hunter**⁽³⁾ - **Julia Jasak**⁽²⁾

SCIEX, Redwood City (1) - SCIEX, Darmstadt, Germany (2) - SCIEX, Concord, Canada (3) - SCIEX, Warrington (4) - Stanford Univ., Palo Alto (5)

Keywords: DIA, SWATH, proteomics, libraries, cloud computing

Introduction

As the use of data independent acquisition (DIA) grows in proteomics research, the need for improved data

processing workflows increases. The most common data processing workflow is to use spectral ion libraries to drive targeted extraction of peptide / fragment areas from the data, using the m/z and retention time information contained in the library. Increasing the size and quality of the ion library has been shown to increase the proteins reliably quantified from a dataset. Retention time (RT) correlation between ion library and the dataset is another key factor that determines quality of data extraction. Simplifying how ion libraries are built and used during data extraction was explored in this work.

Methods

Two algorithms were explored to simplify the data extraction workflow for SWATH® Acquisition data within the OneOmics™ Project cloud processing pipeline, the automatic merging of related ion libraries followed by auto RT calibration. During library merging, ProteinPilot™ Software group files are combined by selecting the largest file as seed library, then peptides from smaller libraries are merged in using a non-linear calibration strategy. New peptides are added to existing proteins and new proteins are added if not present in seed library. During SWATH processing, endogenous peptides are automatically selected across time bins and the best peptide are chosen based on precursor intensity and ID confidence. Best scoring peak groups are used for RT calibration.

Results

Previous data has shown that increasing the size of the ion libraries used provides improved protein quantitation. Moving from ion libraries generated from 1D data dependent acquisitions to libraries generated from 2D fractionation experiments to pan proteome libraries provides a >100% gain in peptides quantified (~60% gain in proteins) to >300% peptide quantified over 1D libraries (>100% gain in protein). Therefore the ability to easily combine and align ion libraries for greater depth of sample interrogation would provide a significant workflow benefit.

To evaluate the quality of the retention time alignment during SWATH processing using the new AutoRT calibration function, it was compared to using the typical RT calibration process using a set of spiked standards peptides. In 33/48 tests, AutoRT approach found similar or more proteins quantified than the standard approach and similar linear fit equations, indicating the new, easier to implement process was as good or better than the currently used standard peptide workflow.

The quality of library merging on SWATH results was explored with a series of small datasets using non-linear retention time alignment. By merging libraries from replicates of a 1D dataset, a 20% gain in peptides was observed. The gain in peptides identified and quantified by SWATH is primarily attributed to an increase in peptide coverage rather than improved ion selection.

Full pipeline was then used on some biological datasets (PBMC and Mouse cell lines) to measure performance, where a series of libraries were available to merge. results will be discussed.

Conclusions

Non-linear retention time alignment was shown to be effective for merging of libraries with identical or differing gradient lengths resulting in better results.

Novel Aspect

Simplifying DIA data processing in the cloud by pipelining ease of use algorithms

MP-172 / THE IMPROVEMENT IN COGNITIVE FUNCTIONS IN 5XFAD ALZHEIMER'S TRANSGENIC MOUSE MODEL VIA NATURAL COMPOUND ADMINISTRATION

Merve Karayel Başar⁽¹⁾ - *İrem Kırış*⁽¹⁾ - *Büşra Gürel*⁽¹⁾ - *Tomasz Mroczek*⁽²⁾ - *Krystyna Skalicka-Woźniak*⁽²⁾ - *Ahmet Tarık Baykal*⁽¹⁾ - *Mehmet Yalçın Günel*⁽³⁾

Acibadem Mehmet Ali Aydınlar University, Medical Biochemistry, Istanbul, Turkey (1) - Medical University of Lublin, Department of Pharmacognosy, Lublin, Poland (2) - Alaaddin Keykubat University, Department of Physiology, Antalya, Turkey (3)

Introduction

Alzheimer's Disease (AD) is a neurodegenerative disorder characterized by memory loss and cognitive decline. (1)The current AD treatment methods aim to prevent cognitive loss. There is no treatment to slow down progression of the disease or stop it. (2)In this study, the effects of new natural compounds on cognitive impairment will be examined and localization of the compounds will be demonstrated by mass-spectrometry based approach.

Methods

The natural compounds were administered at 2 different doses to 12 months old transgenic mice(5XFAD) for 7 days i.p. (n = 8). Non-transgenic littermates (LM) were used as the control group. Morris Water Maze test was performed to measure cognitive impairment. The MALDI-IMS technique was applied to demonstrate whether the natural compounds pass the blood brain barrier and determine the localization of the compound in the brain regions.

Results and Conclusions

Cognitive changes related to learning and memory were measured by Morris Water Maze test after injection of low and high doses of 6 different compounds. When compared to control groups, it was observed that 2 of the 6 compounds positively affected both the learning process and memory. It was demonstrated that the compounds were localized in certain brain regions after they have passed the blood brain barrier by utilizing MALDI-IMS.

Novel Aspect

This is the first study to show the regional localization of the novel native compounds administered on the 5XFAD mouse model.

References

1. Holtzman DM., John CM., Goate A., Alzheimer's Disease: The Challenge of the Second Century, *Apr 6; 3(77): 77sr1* (2011).
2. Van Dam D., De Deyn P. P., Animal models in the drug discovery pipeline for Alzheimer's disease, *Br J Pharmacol*, 164:1285–1300 (2011).

MP-173 / SIMULTANEOUS DUAL PROTEOME ANALYSIS TOWARDS UNDERSTANDING NEUTROPHIL INTERACTION WITH THE HUMAN-PATHOGENIC FUNGUS ASPERGILLUS FUMIGATUS

Thomas Krüger⁽¹⁾ - **Iordana A. Shopova**⁽¹⁾ - **Tobias Müller**⁽²⁾ - **Marcus Dittrich**⁽²⁾ - **Derek J. Mattern**⁽¹⁾ - **Theresia Conrad**⁽³⁾ - **Olaf Kniemeyer**⁽¹⁾ - **Axel A. Brakhage**⁽⁴⁾

Leibniz Institute for Natural Product Research and Infection Biology (Hans Knöll Institute), Department of Molecular and Applied Microbiology, Jena (1) - Julius Maximilians University Würzburg / Biocenter, Department of Bioinformatics and Research Center for Infectious Diseases (ZINF), Würzburg (2) - Leibniz Institute for Natural Product Research and Infection Biology (Hans Knöll Institute), Department of Systems Biology and Bioinformatics, Jena (3) - Leibniz Institute for Natural Product Research and Infection Biology (Hans Knöll Institute), Friedrich Schiller University Jena, Department of Molecular and Applied Microbiology, Jena (4)

Keywords: *Aspergillus fumigatus*, neutrophil extracellular traps, dual proteome, phosphoproteome, host-pathogen interaction

Introduction

Invasive aspergillosis caused by the fungus *Aspergillus fumigatus* is a life-threatening infection in immunocompromised hosts. Neutrophils are indispensable in eliminating the fungal pathogen by various

strategies, i.a. phagocytosis and neutrophil extracellular traps (NETs) [1, 2]. To discover how fungal hyphae activate pathogen-specific responses of neutrophils and to explore fungal counter-defense mechanisms we conducted a quantitative dual proteome study.

Methods

Neutrophils forming NETs against *A. fumigatus* were compared with chemically induced NETs. Hyphae trapped by NETs were compared against those controls grown under equal conditions devoid of neutrophils. LC-MS/MS-based bottom-up proteomics was used to quantify proteins based on the MS2 reporter ion intensities of multiplexed unique peptides labeled with tandem mass tags (TMT). Protein-protein interactions (PPIs) were studied using a network modeling approach.

Results

We were able to quantify 454 differentially abundant proteins of 843 proteins in total on the fungal side and 335 differentially abundant proteins of 1906 proteins in total on the host side. Neutrophils induce fungal proteins involved in the cell wall stress response, toxic secondary metabolite production, and the unfolded protein response. *A. fumigatus* modulates the trace metal ion utilization. Neutrophils and NETs also strongly repress the fungal protein biosynthetic machinery. Analysis of the neutrophil response upon fungal confrontation revealed an altered stoichiometry of alternative and core histones and a modulation of the nuclear envelope. Neutrophils also induce antimicrobial proteins and a repression of proteins involved in calcium-dependent membrane scaffolding. Furthermore, the neutrophil-fungal interaction caused an impairment of the energy metabolism on both sides. PPI network modeling revealed the induction of an autophagy pathway in *A. fumigatus* and a regulatory role of CDK1 during NETosis.

Conclusions

We investigated the responses of human neutrophils to *A. fumigatus* hyphae using a simultaneous dual proteome analysis. The challenging approach provided us with an ample dataset of differentially abundant proteins on both sides. We generated a model of interaction that allowed us to reveal protein interaction partners and regulatory circuits governing the counter response of *A. fumigatus* and the mechanism of hyphal induced NETosis.

Novel Aspect

Our study revealed the fungistatic effects of NETs and the host-pathogen interplay on the molecular level.

References

1. Heinekamp T., Schmidt H., Lapp K., Pahtz V., Shopova I., Koster-Eiserfunke N., Krüger T., Kniemeyer, O., Brakhage A. A., Seminars in Immunopathology, 37:141-152 (2015).
2. Krüger, T., Luo, T., Schmidt, H., Shopova, I., Kniemeyer, O., Proteomes, 3:467-495(2015).

MP-174 / BIOMARKER DISCOVERY OF ATHEROSCLEROTIC STROKE IN INDIVIDUAL SERUM SAMPLE

Jiyeong Lee⁽¹⁾ - **Sora Mun**⁽²⁾ - **Arum Park**⁽²⁾ - **Yoo-Rim Lee**⁽²⁾ - **Ae Eun Seok**⁽²⁾ - **Hyo-Jin Kim**⁽²⁾ - **Yoo-Jin Lee**⁽²⁾ - **Hee-Gyoo Kang**⁽¹⁾

Department of Biomedical Laboratory Science, College of Health Sciences, Eulji University, Seongnam-si, Gyeonggi-do (1) - Department of Senior Healthcare, BK21 Plus Program, Graduate School, Eulji University, Daejeon (2)

Keywords: Biomarker, Proteomics, Stroke, Subtype

Introduction

Ischemic stroke is classified into large artery atherosclerosis (LAA), cardioembolism and small vessel

occlusion. In this study, LAA that is the most common cause of stroke is analyzed compared to healthy subjects. Patients with LAA stroke is approximately 85%. To increase efficacy of diagnosis and recovery, our aim is to discover biomarker of LAA using individual serum of patients.

Methods

We conducted proteomics analysis of serum protein in patients with LAA stroke. 52 patients and 43 healthy subjects are used to identify global proteins and to compare level of serum proteins. Individual samples was analyzed using Sciex 5600 TripleTOF. Identification and quantification of proteins were conducted by information dependent acquisition (IDA) mode and SWATH-MS acquisition.

Results

As a result, differentially expressed proteins by more than 2 fold were 149. Differentially expressed proteins by more than 2-fold was related with GO processes such as regulated exocytosis, exocytosis, secretion by cell, immune response and regulation of response to stress. And, alternative complement pathway were top scored. Among differentially expressed proteins that is related with alternative complement pathway vitronectin, complement factor B, complement component C7, and complement component C5 had the most statistically significance. Vitronectin, complement factor B and complement component C7 were upregulated and complement component C5 was downregulated in patients with LAA stroke.

Conclusions

A proteomics approach discovered associations of LAA stroke with several immune responses compared to healthy subjects. This study highlights four alternative complement pathway-related proteins including vitronectin, complement factor B, complement component C7, and complement component C5. These four proteins contribute to immune responds and B cell and T cell activation and can be potential biomarker to understand LAA stroke.

Novel Aspect

The proteomics analysis of serum of patients with LAA stroke showed alternative 4 immune-related proteins can help to distinguish subtype of stroke from healthy subjects.

References

1. Adams, Jr., B. H. Bendixen., L. J. Kappelle., J. Biller., B. B. Love., D. L. Gordon., E. E. Marsh., Stroke, 24, 35-41 (1993)
2. Rui Liu., Meng-Xian Pan., Jun-Chun Tang., Ya Zhang., Hua-Bao Liao., Yang Zhuang., Dan Zhao., Qi Wan., Neuroimmunol Neuroinflammation, 4, 158-66 (2017)
3. Furie KL., Jayaraman MV., Stroke, 49, 509-510 (2018)

MP-175 / DEVELOPMENT OF A SIMPLE AND EFFICIENT QUANTIFICATION METHOD OF HUMAN FOLLICULAR FLUID (HFF) PROTEOME

Aleksandra Lewandowska⁽¹⁾ - ***Anna Fel***⁽¹⁾ - ***Katarzyna Macur***⁽¹⁾ - ***Paulina Czaplewska***⁽¹⁾ - ***Joanna Liss***⁽²⁾ - ***Krzysztof Łukaszuk***⁽²⁾ - ***Stanisław Ołdziej***⁽¹⁾

Intercollegiate Faculty of Biotechnology UG-MUG, University of Gdansk, Gdansk (1) - INVICTA Fertility Clinic, Gdansk (2)

Keywords: human follicular fluid, proteomics, SWATH-MS, MED-FASP, Total Protein Approach (TPA)

Introduction

Human follicular fluid (hFF) is a potential source of biomarkers of the oocyte quality, which could be assayed in a non-invasive manner prior to oocyte fertilization. In our previous study [1], we investigated the proteome and peptidome of hFF by SWATH-MS methodology [2]. Here, we combine and compare that

methodology with MED-FASP procedure and quantitative MS measurements by TPA [3] to develop an efficient quantification method of hFF proteins.

Methods

Sample of pooled hFFs from patients undergoing the IVF procedure was acquired from the INVICTA Fertility Clinic in Gdansk. Samples after prior ultrafiltration (>10 kDa fraction) or without preliminary treatment were digested by two distinct protocols: i) in-solution trypsin digestion [1]; ii) FASP or MED-FASP [3]. Afterwards, samples were analyzed by two distinct MS methods: i) SWATH-MS [2]; TPA [3].

Results

We have developed a methodology for hFF proteome and peptidome identification and SWATH-MS quantification strategy. To increase the size of the spectral library used for quantification, we applied extensive multistep fractionation consisting of ultrafiltration, optional immunodepletion, and high pH RP-HPLC. We managed to identify 400 proteins with high confidence: 302 in High Molecular Weight Fraction (HMWF, >10 kDa), and 161 proteins in peptidomic Low Molecular Weight Fraction (LMWF, <10 kDa). In order to improve our method, we intend to conduct comparative experiments of two sample preparation procedures: previously used in-solution digestion and FASP or MED-FASP. Our preliminary results suggest that both protocols allow for high confidence discovery of slightly different sets of proteins. Moreover, ultrafiltration prior digestion improved the outcomes. We plan to carry out more experiments to further test both protocols along with two quantification methods: SWATH-MS and TPA.

Conclusions

Our methodology enables a relatively simple quantitative measurement of a high number of proteins in one run. A time and cost efficient fractionation method, such as ultrafiltration significantly improves the qualitative and quantitative analyses of hFF, at the same time allowing for investigations of scarcely studied peptidomic fraction (LMWF). MED-FASP is a promising sample preparation strategy with a high potential of application in hFF proteomics.

Novel Aspect

This is a first work employing MED-FASP and TPA strategy for hFF proteome identification and quantification, as well as their combination with SWATH-MS.

References

1. Lewandowska A. E., Macur K., Czaplewska P., Liss J., Łukaszuk K., Ołdziej S., Journal of Proteomics, In Press (2018), doi: 10.1016/j.jprot.2018.03.010.
2. Gillet L. C., Navarro P., Tate S., Röst H., Selevsek N., Reiter L., Bonner R., Aebersold R., Molecular and Cellular Proteomics, 11(6), O111.016717 (2012), doi: 10.1074/mcp.O111.016717.

MP-176 / A TARGET FISHING AND SHOTGUN PROTEOMICS APPROACHES TO STUDY THE EFFECT OF A NEW DRUG AGAINST M. TUBERCULOSIS

Marcello Manfredi⁽¹⁾ - **Laurent Chiarelli**⁽²⁾ - **Arianna Buzzi**⁽¹⁾ - **Marilina Pasca**⁽²⁾ - **Fabio Gosetti**⁽¹⁾ - **Elettra Barberis**⁽¹⁾ - **Elisa Robotti**⁽¹⁾ - **Emilio Marengo**⁽¹⁾

University of Piemonte Orientale, Department of Sciences and Technological Innovation, Alessandria (1) - University of Pavia, Department of Biology and Biotechnology, Pavia (2)

Keywords: Quantitative proteomics, Target fishing approach, Click chemistry, Drug effect, M. tuberculosis.

Introduction

Mycobacterium tuberculosis (MT) is the causative agent of tuberculosis. The development of more

effective and shorter MT treatments requires the killing of both MT growing and non-growing bacilli. Today, mass spectrometry based proteomics can address challenges that cannot be approached by genomic analysis. The study of the effect of a drug and the molecular mechanisms contributing to disease pathogenesis can be performed using quantitative proteomics.

Methods

A target fishing proteomic approach for the in vitro study of the effect of a new drug, TP053, on MT was performed. The target fishing experiment, which is based on click chemistry, was used to understand if an active metabolite was able to bind a possible target in MT. Moreover, in order to explore the possible pathways affected by the new drug, MT cell lines were treated with TP053 and a shotgun proteomic approach was used to identify modulated proteins.

Results

The TP53 treatment negatively affected the redox reactions, intermediary metabolism and DNA. There is a down regulation of proteins linked to growth processes and generation of metabolites but also a down regulation of carbon and lipid metabolism and ATP generation proteins. All this data are in agreement with the inhibition of aerobic respiration which is a well known NO target.

There is also a modulation of enzymes involved in lipid metabolism and an over expression of many ribosomal proteins, most of them belonging to the 30s ribosome unit.

Using the target fishing approach we were able to identify a direct interaction of TP053 with 30S ribosomal protein S13.

Conclusions

By a shotgun proteomic approach and a target fishing experiment, it has been demonstrated that the activated TP053 is responsible for a global perturbation in the mycobacterial cells through the release of nitric oxide, as well as by a direct interaction with 30S ribosomal protein S13.

Novel Aspect

A target fishing experiment coupled to shotgun proteomic approach allowed to study the effect of a new drug and the interaction drug-proteins.

MP-177 / A COST-EFFECTIVE APPROACH TO PRODUCE ¹⁵N-LABELLED AMINO ACIDS EMPLOYING CHLAMYDOMONAS REINHARDTII CC503

Juan Manuel Marchante-Gayón⁽¹⁾ - *Jesús Nicolás Carcelén*⁽¹⁾ - *Pablo Rodríguez González*⁽¹⁾ - *Luis Valledor González*⁽²⁾ - *María Jesús Cañal Villanueva*⁽²⁾ - *José Ignacio García Alonso*⁽¹⁾

University of Oviedo, Department of Physical and Analytical Chemistry, Oviedo (1) - University of Oviedo, Department of Biology of Organisms and Systems, Oviedo (2)

Keywords: Chlamydomonas reinhardtii, ¹⁵N-labelled amino acids, GC-MS.

Introduction

The synthesis of isotopically labelled compounds from enriched stable isotopes is an expensive and difficult task. Both chemical and biochemical methods have been proposed to produce labelled molecules, but few cost-effective methods have been described [1]. Previously reported biochemical methods to prepare labelled amino acids have only considered their use as hydrolyzed proteins to grow cell cultures or in proteomics [2].

Methods

The aim of this study was to use the microalgae Chlamydomonas reinhardtii to produce, at laboratory scale,

^{15}N -labelled amino acids with a high isotopic enrichment. To do that, a culture medium containing $^{15}\text{NH}_4\text{Cl}$ was used. The labelled proteins biosynthesized by the microorganism were extracted from the biomass. Then, the ^{15}N -labelled amino acids were obtained after a protein hydrolysis with HCl and their isotopic enrichment was measured by GC–MS [3].

Results

No differences in cell growth were observed when they were feeding with natural abundance N or ^{15}N , demonstrating an absence of isotopic effects. On the other hand, the use of the strain CC503 cw92 mt+ and $^{15}\text{NH}_4\text{Cl}$ in the culture medium reduces the incorporation of natural abundance nitrogen, as this strain is not able to incorporate nitrates or other nitrogen containing compounds. Moreover, this wall deficient mutant facilitated the extraction of proteins, leading to higher extraction yields ($32.5\pm 4.0\%$) than those obtained for a wildtype strain 137c ($12.7\pm 3.8\%$). The average isotopic enrichment obtained for the 16 detected proteinogenic amino acids was $99.56\pm 0.05\%$ demonstrating the lack of natural abundance nitrogen incorporation through the whole procedure. The amino acid content ranged from 18–90 $\mu\text{g}/\text{mL}$, depending on the occurrence of the aminoacids present in the proteins. A derivatization procedure based on the use of MTBSTFA was optimized to convert the amino acids into volatile compounds.

Conclusions

Sixteen ^{15}N -labelled amino acids were obtained with an extremely high enrichment from the strain *Chlamydomonas reinhardtii* CC503. As the effective incorporation of ^{15}N has been proved, other labelled bio-molecules should be present in the cells, and this biomass could be employed to isolate ^{15}N -labelled nucleic acids and other metabolites in future experiments. Additionally, multiple labelling could be performed if ^{13}C were supplied to the culture media.

Novel Aspect

An innovative, simple, and fast method to produce isotopically labelled amino acids, which could be used for future quantitative proteomics studies, is presented.

References

1. Lehmann W. D., Mass Spectrometry Reviews, 36, 58-85 (2017).
2. Arsova B. et al., Trends in Plant Science 17(2), 102-112 (2012).

MP-178 / CHARACTERIZATION OF PROTEIN EXTRACTS FROM YEAST CELLS EXPOSED TO FERMENTED SPIRULINA (*ARTHROSPIRA PLATENSIS*) USING HIGH RESOLUTION ORBITRAP MASS SPECTROMETRY

Jasmina Masten^(1,2) - Berta Cillero-Pastor⁽³⁾ - Ronny Mohren⁽³⁾ - Nives Ogrinc^(1,2) - Polona Jamnik⁽⁴⁾

Department of Environmental Sciences, Jožef Stefan Institute, Ljubljana, Slovenia (1) - Jožef Stefan International Postgraduate School, Ljubljana, Slovenia (2) - The Maastricht MultiModal Molecular Imaging Institute (M4I), Maastricht University, Maastricht, Netherlands (3) - Biotechnical faculty, University of Ljubljana, Ljubljana, Slovenia (4)

Keywords: Spirulina, phycobili protein, microbial fermentation, *Lactobacillus plantarum*, Orbitrap

Introduction

Arthrospira is a cyanobacteria with a high nutritional and therapeutic potential. It is a source of essential and non-essential nutrients and bioactive compounds with antiinflammatory and antioxidant activity [1]. Large proportion of *Arthrospira* dry weight comprises phycobili proteins of high bioactive value [2, 3]. Microbial fermentation can contribute to products nutritional and bioactive efficiency due to new metabolite production and transformation [4].

Methods

Arthrospira samples were fermented with *L. plantarum* culture, protein extraction of fermented and non-fermented samples followed. Yeast culture was treated with protein concentrates and later extracted. Protein content and pathway alterations were analyzed by Bradford assay and HF Orbitrap. Liquid chromatography tandem mass spectrometry was used to determine protein changes and Proteome Discoverer for protein identification against the UniProt database.

Results

Exposure of yeast cells to fermented or non-fermented protein extracts affected yeast cell protein content. A higher protein concentration was determined in yeast cells treated with fermented and non-fermented *Arthrospira* extracts compared to control. Furthermore, label free high mass resolution proteomics enabled us to determine protein changes during microbial fermentation in *Arthrospira* water-phase protein extracts.

Conclusions

Arthrospira microbial fermentation using *Lactobacillus plantarum* culture offers great possibilities for protein efficiency and uptake improvement. During lactic acid fermentation proteins are modulated with enhanced digestibility and bioavailability and adding bioactive properties to the food component. This method enables us to develop a food product with an upgraded functional value.

Novel Aspect

Further analysis of lactic acid fermentation of *Arthrospira* is needed to specify its effects on products nutritional and bioactive properties.

References

1. Chopra K., Bishnoi M., Spirulina in human nutrition and health, Gershwin M. E. and Belay A., Eds. Boca Raton: CRC Press, pp. 101–118 (2008).
2. Khan Z., Bhadouria P., Bisen P. S., Curr. Pharm. Biotechnol., vol. 6, pp. 373–379(2005).
3. McCarty M. F., J. Med. Food, vol. 10, pp. 566–570 (2007).
4. Liu J. -G., Hou C. -W., Lee S. -Y., Chuang Y., Lin C. -C., Process Biochem., vol. 46, pp. 1405–1410(2011).

MP-179 / QCONCAT-BASED ABSOLUTE QUANTIFICATION OF THE LYSOSOMAL PROTEOME BY MULTIPLE REACTION MONITORING

***Peter Mosen*⁽¹⁾ - *Roman Sakson*⁽²⁾ - *Thomas Ruppert*⁽²⁾ - *Volkmar Gieselmann*⁽¹⁾ - *Dominic Winter*⁽¹⁾**

Institute of Biochemistry and Molecular Biology, University of Bonn, Bonn, Germany (1) - Zentrum für Molekulare Biologie, University of Heidelberg, Germany (2)

Keywords: QConCat, MRM, Absolute Protein Quantification, Lysosome

Introduction

Absolute quantification of proteins relies on externally quantified peptide standards which are spiked in the sample [1]. For the quantification of large numbers of proteins, the QConCat strategy [2] efficiently provides multiple internal standard peptides concatenated in one proteolytically cleavable artificial fusion protein. We applied this strategy for absolute quantification of the currently known lysosomal proteome using targeted mass spectrometry.

Methods

12 QconCats were expressed in *E. Coli* using stable isotope labeled amino acids and purified via a His tag at their C-Terminus. Tryptic digestion was performed by in gel digestion, in solution digestion in the presence

of 1 M GCL or RapiGest as well as by the FASP [3] and SP3 [4] strategies. The generated peptides were desalted and analyzed using LC-MSMS. Data analysis was performed by Mascot and scheduled MRM transitions were established using Skyline [5].

Results

Based on in house and published large scale proteomics studies, we selected 422 unique tryptic peptides of lysosomal proteins acting as representatives. We were able to cover 143 lysosomal proteins representing 94 % of the known lysosomal proteome [6]. The selected peptides were used to design 12 QconCats which were reverse translated and generated by gene synthesis. After expression in *E. Coli* and His tag purification of the QconCats, we optimized tryptic digestion. We tested in gel digestion, in solution digestion using RapiGestor 1M GCL as well as the FASP and SP3 approaches. RapiGest gave us the lowest number of missed cleavage sites and was therefore used for sample processing. A mixture of all QconCat peptides was analyzed by LC-MSMS using an Orbitrap Velos and QTrap 5500/6500 mass spectrometers. For 417 of the 422 peptides we defined MRM transitions and are currently in the process of assigning specific retention times. In an initial study, we absolutely quantified 48 lysosomal proteins from mouse embryonic fibroblast using the first 4 QconCats.

Conclusions

We successfully designed, expressed and purified 12 QconCats providing internal standard peptides for absolute quantification of the majority of currently known lysosomal proteins. We optimized purification and proteolytic digestion of the QconCats. Assay development for targeted quantification by MRM is almost finalized. Our assay will allow for the first time to perform absolute quantification of the known proteins of a whole organelle.

Novel Aspect

The approach allows for the first time the absolute quantification for almost all proteins of an organelle, the lysosome.

References

1. Bantscheff M., Lemeer S., Savitski MK., Kuster B., *Analytical and Bioanalytical Chemistry*, 404, 939-965 (2012)
2. Pratt JM., Simpson DM., Doherty MK., Rivers J., Gaskell SJ., Beynon RJ., *Nature Protocols*, 1, 1029-1043 (2006)
3. Wiśniewski JR., Zougman A., Nagaraj N., Mann M., *Nature Methods*, 6,359-362(2009)
4. Hughes CS., Foehr S., Garfield DA., Furlong EE., Steinmetz LM., Krijgsveld J., *Molecular Systems Biology*, 10, 1-10 (2014)
5. MacLean B., Tomazela DM, Shulman N., Chambers M., Finney GL., Frewen B., Kern R., Tabb DL., Liebler DC., MacCoss MJ., *Bioinformatics*, 26, 966-968 (2010)
6. Thelen M., Winter D., Braulke T., Giesemann V., *methods in MolecularBiology*, 1594, 1-18 (2017)

MP-180 / VALIDATION STRATEGIES FOR THE PURITY ASSESSMENT OF PEPTIDE STANDARDS BY ISOTOPE DILUTION MASS SPECTROMETRY

Jesús Nicolás Carcelén - Adriana Gonzalez-Gago - Amanda Suarez-Fernandez - Juan Manuel Marchante-Gayón - Pablo Rodriguez-González - José Ignacio García Alonso

University of Oviedo, Analytical and Physical Chemistry, Oviedo

Keywords: Peptide hydrolysis, Method validation, Isotope Dilution, Amino Acid Analysis, Peptide purity assessment.

Introduction

Traceability of protein and peptide quantifications to the International System of units is still a challenge in analytical chemistry. Due to the lack of certified reference materials peptide purity assessment must be carried out through amino acid analysis. This procedure is based on the addition of isotopically labelled aminoacids before the application of a time-consuming acid hydrolysis.

Methods

The only commercially available peptide certified reference materials SRM 998 (Angiotensin I from NIST) and CRM 6901-a (C-peptide from NMIJ) were used for validation purposes. The hydrolysis time was optimized for both materials applying either classical thermal heating or focused microwaves during hydrolysis. Labelled amino acids were spiked to the samples prior hydrolysis with HCl 6N and the peptides were quantified by Isotope Dilution Mass Spectrometry.

Results

The quantification of the peptides was performed by GC-MSMS employing 5 different isotopically labelled amino acids (Valine, Isoleucine, Leucine, Proline and Phenylalanine), which were previously characterized in terms of concentration by reverse isotope dilution using natural abundance certified analogues. The isotopic enrichment of the labelled amino acids was determined as reported elsewhere [1]. The hydrolysis time required to achieve a complete hydrolysis was optimized for SRM 998 and CRM 6901-a when applying thermal heating at 130°C or focused microwaves. Finally, the purity of four isotopically labelled proteotypic peptides of two protein glaucoma biomarkers was assessed applying the optimized conditions. A final validation of the methodology was carried out by quantifying both proteins in human serum samples by isotope dilution and LC-MSMS. An accurate assessment of the peptide purity can be demonstrated if no statistical differences are obtained in the protein concentrations obtained by each labelled peptide.

Conclusions

The application of focused microwaves significantly reduces the time required to achieve a complete hydrolysis of peptide standards. Besides the hydrolysis of peptide certified reference materials, the determination of proteins in real samples by LC-MSMS through proteotypic peptides allow an additional validation of peptide purity assessment within the optimized hydrolysis method.

Novel Aspect

Focused microwaves to decrease peptide hydrolysis time. Validation with reference materials and tryptic peptides applied for absolute quantification of proteins.

References

1. González-Antuña A., Rodríguez-González P., García Alonso J. I., Journal of Mass Spectrometry, 49, 681 (2014).

MP-181 / QUANTITATIVE MICROPROTEOMICS FOR THE CHARACTERIZATION OF CENTRAL AND PERIPHERAL NERVOUS SYSTEM OF THE TWITCHER MOUSE

Davide Pellegrini⁽¹⁾ - ***Ambra Del Grosso***⁽¹⁾ - ***Lucia Angella***⁽²⁾ - ***Nadia Giordano***⁽³⁾ - ***Ilaria Tonazzini***⁽²⁾ - ***Matteo Caleo***⁽³⁾ - ***Marco Cecchini***⁽⁴⁾ - ***Liam McDonnell***⁽⁵⁾

NEST, Scuola Normale Superiore, Pisa (1) - NEST, Nanoscience Institute - CNR, Pisa (2) - Neuroscience Institute, CNR, Pisa (3) - NEST, Scuola Normale Superiore and Nanoscience Institute - CNR, Pisa (4) - Fondazione Pisana per la Scienza ONLUS, Laboratory of Proteomics, San Giuliano Terme (PI) (5)

Keywords: microproteomics, TMT, Krabbe, Twitcher, nervous system,

Introduction

Krabbe disease (KD) is an autosomal-recessive sphingolipidosis caused by a deficiency of galactosylceramide beta-galactosidase, and is characterized by an accumulation of psychosine in the nervous system. In this study, we applied a quantitative microproteomic workflow for the characterization of the central and peripheral nervous system of the Twitcher mouse (a mouse model for KD). The analyses focused on the corpus callosum, motor cortex and sciatic nerves.

Methods

The brains and sciatic nerves of 5 Twitcher and 5 littermate wild type mice were extracted and snap frozen. Cryosectioning and laser capture microdissection was then used to isolate small 1mm² regions of interest from each mouse brain. The workflow combined SP3 digestion, TMT labeling and high-pH fractionation followed by LC-MS/MS analysis using a nanoLC-Orbitrap Fusion system. Data analysis was performed using Proteome Discoverer 2.1, Perseus and WebGestalt.

Results

We performed three TMT-10plex experiments to separately compare the proteomes of the sciatic nerve, corpus callosum and motor cortex of 5 TWI and 5 WT mice. Overall, more than 3000 protein groups were quantified in each of the three datasets. Principal component analysis was performed on the protein expression levels and the TWI and WT mice fell into two distinct clusters, well separated by PC1 in all datasets. A two-sided t-test (FDR=0.01) was performed to compare the expression values of specific proteins. We found 75, 14 and 386 protein groups differently regulated between TWI and WT mice in the corpus callosum, motor cortex and sciatic nerves datasets, respectively. Gene ontology analysis revealed several pathways upregulated in the TWI mice related to inflammatory response (*i.e.* antigen processing and presentation, and complement pathway) and to phagocytosis and lysosomes.

Conclusions

The comparison between TWI and WT mice proteomes highlighted many proteins with distinct expression profiles. Most of the enriched biological processes in the TWI mice were related to immune response, accumulation of lysosomal proteins and reduced nervous system development. These initial results will be further investigated and validated in order to gain new insights into the molecular mechanisms of KD.

Novel Aspect

This work is the first reported mass spectrometry-based in-depth proteomic characterization of the Twitcher mouse proteome to study Krabbe disease.

References

1. de Graaf E. L., Pellegrini D., McDonnell L. A., Journal of Proteome Research, 15, 4722-4730 (2016)

MP-182 / QUANTITATIVE STUDY OF METABOLIC PATHWAYS IN C. ELEGANS USING MS WESTERN

Bharath Kumar Raghuraman - Andrej Shevchenko

Max Planck Institute of Molecular Cell biology and Genetics, Biological Mass spectrometry, Dresden

Keywords: MS Western, C. elegans, metabolic enzymes, systems biology, stoichiometry, isoform interplay

Introduction

Traditionally, metabolic pathways are quantitatively characterized by correlating gene expression profiles to protein abundances. However, owing to poor congruence between them, an alternative approach was needed to resolve this issue.

In *C. elegans*, metabolic transitions are crucial for development and survival. They are tightly regulated by tweaking dynamically the components of the pathways, but very less is known on how absolute protein amounts can contribute to this regulation.

Methods

53 Metabolic enzymes from central carbon metabolism of *C. elegans* were quantified by MS Western method [1]. Unique quantity peptides were selected from sequences of metabolic enzymes and from reference protein BSA and expressed as a single chimera in $\Delta\text{Lys}\Delta\text{Arg}$ double auxotroph *E. coli* strain supplemented with isotopically labelled arginine and lysine. The isotopically labelled chimera, full length BSA & whole worm extracts were co-digested and analyzed by LC-MS. The abundances of tryptic peptides from the labelled chimeric protein and the unlabelled target protein were compared and absolute molar quantities of target proteins were calculated.

Results

53 metabolic enzymes spanning glycolysis, gluconeogenesis, TCA and glyoxylate shunt were accurately quantified in fmole/worm using a single isotopically labeled 300 kDa chimera protein standard comprising, in total 200 quantity peptides. Molar concentrations of endogenous enzymes determined with CV <10% allowed to determine stoichiometric ratios between subunits of multiprotein complexes and between individual members of metabolic pathways. To the best of our knowledge, for the first time this work provided a quantitative overview of the molecular organization of 4 metabolic pathways at the protein level. Although the abundance of enzymes is independently regulated, their molar abundances were ranging from 0.03 fmole/worm to 4 fmole/worm and within each pathway, differed less than 4-fold. High accuracy of the protein quantification allowed to confidently determine concentration difference of 2-fold, typically disregarded by label free quantification. By comparing molar abundances of metabolic enzymes in several mutants, this work provided molecular insight into the metabolic regulation at different stages of *C. elegans* development, including the dormant dauer stage.

Conclusions

Absolute quantities of 53 metabolic enzymes provided crucial evidences of tight regulation of metabolic pathways. The changes in the specific gluconeogenic enzymes primed towards the conclusion that dormant worm uses gluconeogenesis as their survival strategy by exploiting all possible metabolic bypasses [2].

Novel Aspect

Comprehensive absolute quantitation of 53 metabolic enzymes spanning 4 major pathways in one experiment with a CV < 10 %. Detailed study of isoform specificity in pathway regulation. Easy estimation of subunit stoichiometry of enzyme complexes. Unique repository for system biologists for pathway studies.

References

1. Kumar. M et al., Molecular & Cellular Proteomics, 17, 384-396 (2018).
2. Erkut Cihan, Vamshidhar R Gade, Sunil Laxman, and Teymuras V Kurzchalia., eLife 5: 1–24 (2016)

MP-183 / APPLICATION OF MANUAL CHROMATOGRAPHY DEVICE FOR FRACTIONATION OF COMPLEX PEPTIDE SAMPLES PRIOR NANOESI-MS/MS

Pavel Rehulka⁽¹⁾ - **Helena Rehulkova**⁽¹⁾ - **Vojtech Franc**⁽²⁾ - **Marek Sebela**⁽²⁾ - **Rudolf Kupcik**⁽³⁾ - **Zuzana Bilkova**⁽³⁾ - **Jiri Stulik**⁽¹⁾

Faculty of Military Health Sciences, University of Defence, Hradec Kralove, Czech Republic (1) - Faculty of Science, Palacký University, Olomouc (2) - Faculty of Chemical Technology, University of Pardubice, Pardubice (3)

Keywords: microgradient device, peptide sample fractionation, basic reversed phase chromatography, nanoLC-MS/MS

Introduction

Proteomic analyses using mass spectrometry often rely on suitable chromatography separation of complex peptide mixtures. In some cases, multidimensional chromatography application is even needed for a deeper insight of proteome. A combination of basic reversed phase chromatography fractionation with subsequent acidic reversed phase chromatography coupled to ESI-MS/MS analysis may be applied.

Methods

Manual gradient chromatography device was used for fractionation of complex peptide mixture. Mobile phases were composed of acetonitrile in 20 mM ammonium formate at pH 10. Microcolumn (250 µm i.d., 50 mm length) contained core-shell particles (Kinetex EVO C18, 2.6 µm) and sample was fractionated either into 17 or 33 fractions that were further dried and analyzed using nanoLC-MS/MS (Q-Exactive).

Results

The nanoLC-MS/MS analysis of a complex peptide mixture without fractionation enabled identification of 1975 protein groups. The fractionation of this sample using manual chromatography device in basic reversed phase chromatography revealed identification of up to 6197 protein groups. Not only the number of identified proteins was increased, but also the number of identified peptide sequences was increased about 4 times.

Conclusions

The manual chromatography device fractionates low amounts of peptides (tens of micrograms) into several fractions that provide a much deeper insight into the analyzed proteome after the nanoLC-MS/MS analysis. The device is advantageous for its simplicity of application, speed of fractionation (about 45 min per sample) and its low cost.

Novel Aspect

Manual chromatography device for peptide sample fractionation prior to nanoLC-MS/MS analysis.

MP-184 / DISCUSSION OF THE PROS AND CONS OF ISOBARIC LABELLING COMPARED TO SINGLE SHOT DATA-INDEPENDENT ACQUISITION

Jan Muntel⁽¹⁾ - *Roland Bruderer*⁽¹⁾ - *Joanna M Kirkpatrick*⁽²⁾ - *Oliver M. Bernhardt*⁽¹⁾ - *Lynn Verbeke*⁽¹⁾ - *Tejas Gandhi*⁽¹⁾ - *Ting Huang*⁽³⁾ - *Olga Vitek*⁽³⁾ - *Alessandro Ori*⁽²⁾ - *Lukas Reiter*⁽¹⁾

Biognosys, R&D, Schlieren (1) - Leibniz Institute on Aging, Fritz Lipmann Institute, Jena (2) - College of Computer and Information Science, Northeastern University, Boston (3)

Keywords: proteomics, Isobaric labelling, data independent acquisition, quantification, comparison

Introduction

Label free quantification (LFQ) and isobaric labelling quantification (ILQ) are the two most popular quantification workflows in discovery proteomics. Here, we wanted to objectively compare the two workflows given a fixed amount of LC-MS time. For ILQ, we chose TMT 10-plex, for LFQ we chose single shot data-independent acquisition (DIA). We joined forces with a laboratory at FLI where ILQ is well established to compare the two workflows.

Methods

We generated a gold standard set of ten samples composed of a complex mouse cerebellum background spiked with the UPS2 protein standard (Sigma) in five different concentrations. LC-MS data was acquired on

an Orbitrap Fusion Lumos (Thermo) in both laboratories and with both workflows. For TMT an MS3 method was used, for DIA a method as described in (Bruderer et al. 2017) was chosen. Data were analyzed with Proteome Discoverer (Thermo) and Spectronaut (Biognosys).

Results

In the four experiments (two workflows in two laboratories) we identified 48,500 (FLI TMT), 47,700 (BGS TMT), 51,800 (FLI DIA) and 53,500 (BGS DIA) peptides. On the protein level this corresponded to 6,500 (FLI TMT), 6,200 (BGS TMT), 5,700 (FLI DIA) and 5,900 (BGS DIA) identifications. 16 (FLI TMT), 15 (BGS TMT), 14 (FLI DIA) and 17 (BGS DIA) UPS2 spike in proteins were detected.

Generally, protein identifications were higher with TMT and protein sequence coverage higher with DIA. The quantitative performance in detection of differential abundance was very similar with a slight edge towards LFQ-DIA at both sites (FLI and BGS) and for two criteria used (partial area under the curve or true positives in the top 200 candidates).

We also performed an ILQ-TMT experiment with very deep pre-fractionation (24 fractions). This illustrated that ILQ-TMT is particularly suited for very deep proteome coverage. In this experiment we were able to quantify 8,800 protein groups, 160,000 precursors and 26 of the UPS2 proteins.

Conclusions

The quantitative performance of both methods was very similar in both laboratories. It turned out that for the analysis of the DIA data it was beneficial to use both quantitative information levels available in DIA, at MS1 and MS2 level. We see the strength in LFQ DIA particularly in a single shot workflow for dozens to thousands of samples and for ILQ TMT particularly for very deep coverage of a few samples.

Novel Aspect

Comparison of isobaric labelling workflow with single shot label free DIA quantification on a gold standard sample set

MP-185 / HIGH SPEED, HIGH SENSITIVITY AND HIGHLY REPRODUCIBLE AND ACCURATE LABEL FREE QUANTIFICATION USING THE PASEF METHOD ON A TIMS QTOF

Gary Kruppa⁽¹⁾ - **Markus Lubeck**⁽²⁾ - **Heiner Koch**⁽²⁾ - **Verena Tellström**⁽³⁾ - **Paul Shan**⁽⁴⁾ - **Jürgen Cox**⁽⁵⁾ - **Scarlet Koch**⁽²⁾

Bruker Daltonics Inc, Proteomics, Billerica (1) - Bruker Daltonik GmbH, Proteomics, Bremen (2) - Bruker Daltonik GmbH, Application, Bremen (3) - Bioinformatics Solutions inc, Development, Waterloo (4) - Max Planck Institute für Biochemie, Computational Systems Biochemistry, Martinsried (5)

Keywords: TIMS, Targeted Proteomics, PASEF, DIA

Introduction

Mass spectrometry-based proteomics has become a powerful technology for the identification and quantification of proteins. We recently introduced the Parallel Accumulation Serial Fragmentation (PASEF) method on TIMS QTOF instrument, delivering high robustness, speed and sensitivity, for shotgun proteomics. Here we show reproducible and accurate quantification using the PASEF acquisition method and the optimized PEAKS and MaxQuant software packages.

Methods

The peptide mixtures (< 200 ng) were separated with a nanoElute (Bruker Daltonics) nano-flow LC on-line coupled to a high-resolution TIMS QTOF (Bruker Daltonics). Chromatographic separation was carried out at a flow rate of 400 nl/min using a 30-90 min linear gradient. LC-MS/MS data were acquired using the PASEF method. Data analysis was performed using PEAKS studio (Bioinformatics Solutions Inc.) and MaxQuant (Jürgen Cox, Max Planck Institute of Biochemistry).

Results

200 ng of a complex peptide mixture derived from a mammalian cell line in the single-run format were analyzed using a 90 min gradient and optimized PASEF parameters. More than 5300 proteins families were identified in each run with a good reproducibility when comparing the protein identities between replicates. Comparing label free intensities between technical replicates shows excellent reproducibility with a $R^2 = 0.98$ using either a 90 min or with a 30 min gradient. To evaluate the reproducibility of PASEF, we determined the coefficient of variation of the label free intensities. Also here good reproducibility could be shown. For the evaluation of the accuracy of quantification over a wide concentration range we performed a spike-in experiment. We spiked E. coli digest in different ratios (1:2, 1:5 and 1:10) into HeLa digest and injected 100 ng sample material and measured each concentration in triplicate analyses. With all different concentrations, expected ratios could be determined.

Conclusions

tims-PASEF allows to perform reproducible, accurate and high sensitive label-free quantitation from shotgun proteomics samples.

Novel Aspect

Label free quantification using PASEF

MP-186 / FAIMS REDUCES INTERFERENCE AND IMPROVES QUANTITATIVE ACCURACY IN MULTIPLEXED PROTEOMICS

Devin Schweppe⁽¹⁾ - *Satendra Prasad*⁽²⁾ - *Romain Huguet*⁽²⁾ - *Jose Navarrete Perea*⁽³⁾ - *João Paulo*⁽³⁾ - *Mark Jedrychowski*⁽³⁾ - *Graeme McAlister*⁽²⁾ - *Susan Abbatiello*⁽²⁾ - *Michael Belford*⁽²⁾ - *Vlad Zabrouskov*⁽²⁾ - *Steven Gygi*⁽³⁾

Harvard Medical School, 240 Longwood Ave, Building C, Boston (1) - Thermo Fisher Scientific, Life Sciences Mass Spectrometry, San Jose (2) - Harvard Medical School, Cell Biology, Boston (3)

Keywords: Multiplexed, TMT, quantitative proteomics, FAIMS

Introduction

Ratio compression reduces dynamic range in multiplexed quantitative proteomics studies, particularly for high-resolution MS2 methods. SPS-MS3 has been highly successful for improving quantitative accuracy by reducing interference-based compression. Building on this, we tested whether the use of high-Frequency Asymmetric Ion Mobility Spectrometry (FAIMS) could further reduce interference in multiplexed quantitative analyses on an Orbitrap Fusion Lumos.

Methods

Three complex samples were used to test the utility of FAIMS to reduce interference: (1) Triple-knockout yeast (TKO, 100% yeast), (2) canonical interference standard (CIS, 46% yeast), and (3) a novel small change interference model (SCI, 10% yeast). Known peptide ratios across channels ranged from 1.5-fold to infinite. Peptides were analyzed on an Orbitrap Fusion Lumos (Thermo Scientific) for analysis with or without a FAIMS source (three CVs tested: -40, -60, -80V).

Results

Thermo recently developed a FAIMS source with improved measurement sensitivity/speed and FAIMS was previously shown to reduce interfering precursors by ion-mobility-based separation [1]. Here, we tested if Thermo's new FAIMS source could reduce co-eluting interfering ions to improve quantitative accuracy for multiplexed samples.

We compared HRMS2 with and without FAIMS using a CIS (labeled yeast dilutions with/without mouse background), and strikingly observed a 55% reduction in quantitative interference. We identified 2288 total yeast peptides with FAIMS-HRMS2 compared to 2423 total yeast peptides for HRMS2 without FAIMS (94%) indicating similar PSM sensitivity with or without FAIMS.

For the TKO sample [2], a 9-plex with triplicates of three KO strains (Δ met6, Δ ura2, Δ pkf2), FAIMS improved the average interference free index by 22% (reduced avg. interference by 62%) compared to the same method without FAIMS. Additionally, we observed a 2.2-fold increase in total peptide identifications for FAIMS-HRMS2 compared to FAIMS-SPS-MS3.

Conclusions

We observed that FAIMS-SPS-MS3 (with the SCI model) had the least quantitative interference (most accurate ratios), followed by non-FAIMS-SPS-MS3 and closely thereafter by FAIMS-HRMS2. Non-FAIMS-HRMS2 had the highest levels of interference and ratio compression. Total peptide identifications for FAIMS-HRMS2 were 54% higher than for SPS-MS3. Our findings suggest that Thermo's new FAIMS source would be a valuable tool for any multiplexed quantitative workflow.

Novel Aspect

FAIMS (HRMS2 or SPS-MS3) consistently improved quantitative accuracy without sacrificing peptide identifications for multiplexed samples compared to standard methods.

References

1. Pfammatter S., Bonneil E., Thibault P., Journal of Proteome Research, 15, 4653-4665 (2016).
2. Paulo J. A., O'Connell J. D., Gygi S. P., JASMS, 27, 1620-1625 (2016).

MP-187 / A TARGETED MULTIPLE REACTION MONITORING METHOD FOR QUANTIFYING NAV1.6 PROTEIN IN HEK293 CELLS AND MOUSE HIPPOCAMPUS

M. Tung^(1,2) - **R. Kwan**⁽²⁾ - **C. Xi**⁽³⁾ - **J. Li**⁽³⁾ - **C. Dang**^(1,3) - **L. Sojo**^(1,2)

Compound Properties Group, Bioanalytical Laboratory, Xenon Pharmaceuticals Inc., Burnaby, B. C., Canada (1) - Department of Chemistry, Simon Fraser University, Burnaby, B. C., Canada. (2) Department of Chemistry University of Victoria, Victoria, B. C., Canada (3)

Keywords: Nav1.6, targeted proteomics, MRM, Human Embryonic Kidney (HEK) 293, Scn8a hippocampus

Introduction

Nav1.6 expression and activity in hippocampus are shown to mediate seizure resistance (1). HEK293 cells expressing human sodium channel Nav1.6 are used to screen compounds targeting Nav1.6 for the discovery of drug candidates to treat epilepsy. Expression of Nav1.6 proteins is traditionally analyzed by Western blot analysis. Here we show that targeted multiple reaction monitoring mass spectroscopy can provide protein quantitation using stable isotope dilution of unique Nav1.6 peptides.

Methods

Membrane proteins were extracted from Nav1.6 expressing HEK293 cells, and from Scn8aN1768D/+ and WT mouse hippocampus (2) using Mem-PER Plus kit (Thermo 89842). Detergents were removed using spin columns (Pierce 87777). Protein (80 μ g) from Nav1.6 HEK293 cells and 240 μ g of Scn8a and WT hippocampal protein were subjected to in-solution tryptic digestion (Thermo 89895). Samples were dried and reconstituted in 50 μ L 1:1 ACN:water + 0.1% TFA for UHPLC-MS/MS analysis.

Results

A total of 39 tryptic peptides, based on an in-silico digestion of a complete human Nav1.6 protein were targeted for MRM analysis of Nav1.6 expressed in a stably transfected HEK293 cells as well as in mouse hippocampus. Out of the 39 peptides, 11 were detected with at least three specific fragment ions. The remaining peptides were present, but not all fragments were detected. From these 11 peptides, seven were unique to Nav1.6. The other four peptides shared homology with other sodium channels. Peptide DSLFIPIR (480.781 m/z , 598-605) was selected for Nav1.6 quantification due to its intense and robust signal. Light and heavy (internal standard) versions of DSLFIPIR were synthesized and used to quantify Nav1.6 using an isotopic dilution method. Nav1.6 HEK293 cells showed Nav1.6 expression of 3.3 fmol/ μg of total membrane protein. Scn8a^{N1768D/+} and WT mouse hippocampus showed 0.41 and 0.39 fmol/ μg of total membrane protein, respectively.

Conclusions

A simple and fast turnaround proteomic work-up was tested for detecting and quantifying Nav1.6 expression in a stably transfected Nav1.6 HEK293 cell line as well as in Scn8a mouse hippocampus. The success of these Nav1.6 proteomic studies enables the use of proteomics to screen and select optimal cell lines for epilepsy drug development. Changes in Nav1.6 expression in mouse hippocampus could also be monitored using proteomics. The same approach can be applied to other ion channels of interest.

Novel Aspect

This is the first time that quantification of Nav1.6 protein expression in HEK293 cell lines and mouse hippocampus is reported.

References

1. Christopher D Makinson et al., Role of the hippocampus in Nav1.6 (Scn8a) mediated seizure resistance, *Neurobiol. Dis.* 2014, 68, 16–25.
2. Jacy L. Wagnon et al., Convulsive seizures and SUDEP in a mouse model of SCN8A epileptic encephalopathy, *Hum. Molec. Gen.* 2015, 24, 506–515.

MP-188 / COMPARATIVE PROTEOGENOMIC ANALYSIS OF COLISTIN-DEPENDENCE IN ACINETOBACTER BAUMANNII

Sherley Chamoun - Jenny Welander - Carina Claesson - Elena Vikström - Maria Turkina

Linköping University, Department of Clinical and Experimental Medicine, Linköping

Keywords: proteogenomics, label free quantification, multidrug-resistant bacteria, hospital-acquired infections

Introduction

The multi-drug resistant bacteria *A. baumannii* belong to the group of human opportunistic pathogens responsible for most common life threatening nosocomial infections. The resistance of *A. baumannii* to the last-resort antibiotic colistin is increasing and, recently, a high rate of colistin dependence has been demonstrated [1, 2]. However, the mechanisms of antibiotic resistance and dependence are not yet fully understood.

Methods

Colistin-susceptible clinical isolate of *A. baumannii* and its subpopulation that exhibited a colistin-dependent growth pattern were studied. The whole-genome sequencing (WGS) and the shotgun proteome analysis by high-resolution mass spectrometry with following bioinformatics analyses were performed to characterize the molecular mechanisms of antibiotic dependence in *A. baumannii*.

Results

We utilized a proteogenomic approach to explore the molecular mechanisms of colistin dependence. The bioinformatics analysis of the obtained WGS data identified multiple antimicrobial resistance genes which were identical in colistin-susceptible and colistin-dependent subpopulations; however colistin-resistant genes were not detected by ResFinder server. WGS data were used to construct a custom strain-specific six-reading frame database, which was used for the proteomic analysis. We found that 183 proteins were expressed differently in the two phenotypes and 119 of them were downregulated in the colistin-dependent subpopulation. The significant part of these proteins are metabolic enzymes, chaperones, outer membrane proteins, efflux pumps and porins. Notably, in addition to novel proteins, several of differentially expressed proteins are known to be involved in colistin-resistance mechanisms.

Conclusions

Our study shows that the significant part of proteome is modified in colistin-dependent sub population of *A. baumannii*. This proteogenomic study also indicates that the colistin-dependence mechanisms are closely related to colistin-resistance. Our research provides a new analytical platform and a foundation for further research of multiresistant *A. baumannii*.

Novel Aspect

This is the first study to use proteogenomics approach for identification of the molecular mechanisms behind antibiotic dependence in *A. baumannii*.

References

1. Nordqvist H., Nilsson LE., Claesson C., Eur J Clin Microbiol Infect Dis. 35(11),1845-1850, (2016)
2. Lee JY., Chung ES., Ko KS., Sci Rep. 7(1), 14216, (2017).

MP-189 / TMT QUANTIFICATION USING SPECTRAL LIBRARIES IN SPECTROMINE

Lynn Verbeke⁽¹⁾ - **Tejas Gandhi**⁽¹⁾ - **Oliver M. Bernhardt**⁽¹⁾ - **Lukas Reiter**⁽²⁾

Biognosys AG, Bioinformatics, Schlieren (1) - Biognosys AG, R&D, Schlieren (2)

Keywords: TMT-quantification, spectral library searching, SpectroMine

Introduction

SpectroMine can search both data-dependent (DDA) and data-independent acquisition (DIA) data using in silico generated search spaces (from FASTA files) or spectral libraries. The advantage of using libraries is the greatly reduced search space, leading to a speed-up and improved sensitivity/reproducibility [1]. Here we set out to see if a library-based targeted approach can lead to improvements, both in identification and quantification, in a TMT-experiment.

Methods

We used a published data set of 7 TMT 10-plex experiments, all from different mouse tissues. Every experiment consists of 12 fractions of 5 controls and 5 nicotine-treated samples [2]. We combined all 7 experiments to generate a library in SpectroMine with 1% false discovery rate at PSM, peptide, protein group level. SpectroMine also identifies the TMT reporter ions, in MS2 or MS3 scans. These are used to calculate differential abundances between conditions.

Results

The generated library has a size of 673'204 fragments, 112'587 precursors, 10'120 protein groups (11 hours for processing 84 runs). We used this library to search the brain experiment. In the case of library searching, SpectroMine uses normalized retention time (iRT) [3] available from the libraries to further improve the search: only peptides within an expected retention time window will be matched to the scans. This resulted

in 48'126 precursors and 7'627 protein groups (20 min). With a Q-value of maximum 0.05 and a minimal absolute log ratio, we found 18 candidates. When searching this data with the mouse FASTA file, we find 39'098 precursors, 6'617 protein groups (2 hours). Applying the same criteria, leads to only 16 candidates. Comparing the quantities of the candidates we obtained with those of the paper, resulted in correlations of 0.98 for both the FASTA as well as for the library approach.

Conclusions

When using a library instead of a FASTA file to identify TMT-labeled peptides, we managed to get an increase of 23% at precursor level and of 15% at protein group level. On top of that, quantification can still be applied by extracting the reporter ions from the MS3 scans. Very similar quantification results are obtained, showing that the library-searching approach of DDA-data can be combined with TMT-quantification.

Novel Aspect

Library searching of TMT-labeled data, while making use of normalized retention time (iRT) to increase identifications while maintaining quantitative accuracy.

References

1. Lam H, Deutsch EW, Eddes JS, Eng JK, King N, Stein SE, Aebersold R, Proteomics, 7, 655 (2007).
2. Joao A. Paulo, Mark P. Jedrychowski, Edward T. Chouchani, Lawrence Kazak, Steven P. Gygi, Proteomics, 10, 1002 (2018).
3. Roland Bruderer, Oliver M. Bernhardt, Tejas Gandhi, Lukas Reiter, Proteomics, 16, 2246 (2016).

MP-190 / MULTIPLEX QUANTITATIVE ASSAY OF INFLAMMATION PROTEIN MARKERS IN FECES SAMPLES

Veronika Vidova⁽¹⁾ - Vojtech Thon⁽²⁾ - Jana Klánová⁽¹⁾ - Zdenek Spacil⁽¹⁾

Reaearch Centre (RECETOX), Masaryk University, Faculty of Science, Brno (1) - St. Anne's University Hospital, Department of immunology, Brno (2)

Keywords: inflammation markers, feces, selected reaction monitoring, targeted proteomics

Introduction

Inflammatory bowel disease (IBD) and food allergies represent a wide group of intestinal disorders that may lead to chronic inflammation of gastrointestinal tract (*e.g.* Crohn's disease, celiac disease). Food allergies are abnormal reactions to an allergen presented in food however the cause of IBD remains unclear [1,2]. Quantification of specific inflammation markers in feces allow to characterize an inflammation and estimate its progress in digestive tract.

Methods

Proteins were extracted from both lyophilized and fresh feces samples for comparison, using 1x phosphate saline buffer and centrifugation. The supernatant was subjected to bottom up proteomics workflow. We have tested isolation of neutrophils and eosinophils from feces and characterize their CD markers. Proteotypic tryptic peptides were selected, inspected for isoforms and SNPs and quantified using UHPLC-SRM analysis.

Results

We have developed and validated multiplex protein assay based on ultra-high performance liquid chromatography (UHPLC) and selected reaction monitoring (SRM) mass spectrometry for the absolute quantification of six inflammatory proteins (*i.e.* ECP, EDN, CAL S100A8, CAL S100A9, F-MPO and A1AT) in feces. Inflammation proteins were quantified in both lyophilized and fresh feces samples. Higher

concentrations of inflammation proteins were quantified in lyophilized samples, most likely due contribution from lysis of neutrophils and eosinophils, which contain these inflammation proteins and occur in feces. Lyophilization allowed us to normalize concentration of each protein to dry weight of lyophilized feces and to achieve lower limit of quantification. Certain inflammation proteins were cross-quantified using ELISA and results of both methods were correlated.

Conclusions

We have developed UHPLC-SRM method for quantitative analysis of selected inflammation proteins in feces to identify an inflammatory status and potentially a specific type of inflammation. We have cross-validate quantitative results of newly developed UHPLC-SRM assay to the results obtained using ELISA assay as a reference method.

Novel Aspect

A novel multiplex method applicable to study development of gastrointestinal immune system, development of IBD disease and food allergies.

References

1. Lopez, R. N. ; Leach, S. T. ; Lemberg, D. A. ; Duvoisin, G. ; Gearry, R. B. ; Day, A. S. J. Gastroenterol. Hepatol., 32, 577–582 (2017).
2. Ohtsuka, Y. Pediatr. Int., 57, 22–29 (2015).
3. Magnusson, J. ; Gellerstedt, M. ; Ahlstedt, S. ; Andersson, B. ; Bengtsson, U. ; Telemo, E. ; Hansson, T. ; Peterson, C. G. B. Clin. Exp. Allergy, 33, 1052–1059 (2003).

MP-191 / DEVELOPMENT OF A QUALITY CONTROL STANDARD FOR TANDEM MASS TAGS (TMT) WORKFLOWS

Jae Choi⁽¹⁾ - *Aaron Robitaille*⁽²⁾ - *Tabiwang N. Arrey*⁽³⁾ - *Rosa Viner*⁽²⁾ - *Andreas Huhmer*⁽²⁾ - *John C. Rogers*⁽¹⁾

Thermo Fisher Scientific, research and development, Rockford (1) - Thermo Fisher Scientific, marketing, San Jose (2) - Thermo Fisher Scientific, marketing, Bremen (3)

Keywords: TMT, proteomic standard, QC, method development

Introduction

Quantitative proteomics strategies using Tandem Mass Tags (TMT) enable sample multiplexing and precise measurement of protein abundance. However, co-isolated ion interference can suppress accurate ratio quantification and thereby mask true differences in protein abundance across biological systems [1]. To be able to detect and diagnosis co-isolation interference, we developed the Thermo Scientific™ TMT11plex™ labeled peptide reference sample.

Methods

Lysates from four strains of *Saccharomyces cerevisiae* were prepared and labeled with TMT11plex reagents and mixed in equimolar ratios. Peptides from each knockout strain were labeled in triplicate, while the parental line was labeled in duplicate. The samples were then analyzed on a range of Thermo Scientific™ Orbitrap™ instruments using LC-HCD MS2 or MS3 fragmentations. Data analysis was performed with Thermo Scientific™ Proteome Discoverer™ 2.2 software.

Results

To develop a widely applicable TMT11plex standard, we modified the selection of yeast strains from Paulo et al. 2016 [1] and included the parental strain BY4742 as reference channels. We show the yeast triple knockout (TKO) standard has less than 20% variability with lot to lot reproducibility. We then used the

TMT11plex standard to establish a standardized workflow for liquid chromatography, MS acquisition on Orbitrap mass spectrometers, and data analysis software. Synchronous precursor selection (SPS) based methods provided the best accuracy and precision as compared to MS2 methods. Depending on precursor isolation purity and notches (fragments) selected for MS3, ratios for knock out proteins can be still distorted. We then implemented an additional filter into the Proteome Discoverer software quantitation node to exclude precursors if more than the user defined number of isolated fragments were not from the identified peptide. Implementation of a new filter significantly improved quantitation accuracy and precision for all proteins.

Conclusions

We demonstrate that the TMT11plex standard can be used as an LC-MS system suitability standard to measure and optimize protein/peptides identification, acquisition, and data analysis methods to limit co-isolation interference, while also diagnosing MS instrument status by monitoring mass accuracy, ion injection time, and reporter ion signal to noise. The TMT11plex TKO standard provides users a tool to measure and assess different mass spectrometry approaches.

Novel Aspect

We present a widely-applicable, high quality TMT11plex yeast peptide standard for LC-MS quality control and application development.

References

1. Paulo J. A., O'Connell J. D., Gygi S. P. *J Am Soc Mass Spectrom.* 27(10): 1620–1625, (2016)

MP-192 / SOLVING THE YELLOW MYSTERY OF PAPAVER NUDICAULE WITH AN INTEGRATED -OMICS APPROACH

Bettina Dudek⁽¹⁾ - *Natalie Wielsch*⁽²⁾ - *Riya Christina Menezes*⁽²⁾ - *Anne-Christin Warskulat*⁽¹⁾ - *Yvonne Hupfer*⁽²⁾ - *Sybille Lorenz*⁽²⁾ - *Heiko Vogel*⁽³⁾ - *Christian Paetz*⁽¹⁾ - *Bernd Schneider*⁽¹⁾ - *Aleš Svatoš*⁽²⁾

Max Planck Institute for Chemical Ecology, NMR Group, Jena (1) - Max Planck Institute for Chemical Ecology, Mass Spectrometry / Proteomics Group, Jena (2) - Max Planck Institute for Chemical Ecology, Department of Entomology, Jena (3)

Keywords: biosynthesis, integrated -omics approach, nudicaulins, Papaver nudicaule

Introduction

Indole alkaloids are a widespread group of secondary metabolites with huge structure variety. Nudicaulins, yellow flower pigments from Papaver nudicaule (Island poppy), are unique representatives of this compound class, combining structural elements from indole and flavonoids. [1] Previous studies identified indole and pelargonidin glucosides as final precursors. [2] The aim of our work is to understand the biosynthesis of the nudicaulins.

Methods

Our approach on the biosynthesis of nudicaulins is based on four pillars: NMR studies on their chemistry, UPLC-HRMS and LC-UV/ Vis based metabolomics, transcriptomics, and DIGE-based quantitative proteomics.

Results

The development of the P. nudicaule flowers was divided in five stages. These stages were monitored by transcriptomics, proteomics, and metabolomics approaches. Thereby, the most important genes and their products involved in the phenyl propanoid/polyketide as well as indole biosynthetic pathway were identified. Relevant metabolites were semi-quantified over the flower development.

Additional experiments showed that the final nudicaulin formation from indole and pelargonidin glucosides is also possible in vitro under acidic pH conditions. This raises the question about the exact mechanism of this reaction in vivo, which is a work in progress.

Conclusions

Based on transcriptomic, proteomic, metabolomic, and NMR studies we propose a biosynthetic pathway of nudicaulins in yellow petals of *P. nudicaule*. This opens new insights into this unique class of indole alkaloids.

Novel Aspect

With the help of a multi method approach, our investigation provides access to the unique biosynthetic pathway of nudicaulins, the yellow pigments of *P. nudicaule* petals.

References

1. Tatsis E. C. ; Schaumlöffel A. ; Warskulat A. -C. ; Massiot G. ; Schneider B. ; Bringmann G. *Org. Lett.*, 15, 156–159 (2013).
2. Warskulat A. -C. ; Tatsis E. C. ; Dudek B. ; Kai M. ; Lorenz S. ; Schneider B. *ChemBioChem*, 17, 318–327 (2016).

MP-193 / A STREAMLINED WORKFLOW FOR HIGH-THROUGHPUT, PRECISE, AND COMPREHENSIVE LARGE-SCALE QUANTITATIVE PROTEOMICS ANALYSIS

Yue Xuan⁽¹⁾ - **Yue Zhou**⁽²⁾ - **Sebastien Gallien**⁽³⁾ - **Oleksandr Boychenko**⁽⁴⁾ - **Joshua J. Nicklay**⁽⁵⁾ - **Jenny Ho**⁽⁶⁾ - **Scott Peterman**⁽⁷⁾ - **Pedro Navarro**⁽¹⁾ - **Emily Chen**⁽⁷⁾ - **Ken Miller**⁽¹⁾

Thermo Fisher Scientific, Chromatography and Mass Spectrometry Division (CMD), Bremen (1) - Thermo Fisher Scientific, Chromatography and Mass Spectrometry Division (CMD), Shanghai (2) - Thermo Fisher Scientific, Chromatography and Mass Spectrometry Division (CMD), Paris (3) - Thermo Fisher Scientific, Chromatography and Mass Spectrometry Division (CMD), Germering (4) - Thermo Fisher Scientific, Chromatography and Mass Spectrometry Division (CMD), New Jersey (5) - Thermo Fisher Scientific, Chromatography and Mass Spectrometry Division (CMD), Hemel Hempstead (6) - Thermo Fisher Scientific, Chromatography and Mass Spectrometry Division (CMD), Cambridge (7)

Keywords: label-free quantitation, data-independent acquisition (DIA), high-throughput, reproducibility, robustness

Introduction

A reproducible and robust global proteome profiling workflow is required to confidently determine putative protein marker panels. In this work, we illuminate the impact of key components of the DIA workflow, including chromatography, MS resolution & accuracy, spectral library quality, as well as the statistical analysis to demonstrate the label-free quantitation advantages of the HRMS1 – DIA workflow on a complex multi-proteome sample set.

Methods

Experiments were performed on a Thermo Scientific™ Q Exactive™ HF mass spectrometer using a method consisting of full scan MS with interspersed DIA acquired for qualitative analysis. Chromatographic separations were performed using a Thermo Scientific™ UltiMate™ 3000 RSLC system, equipped for capillary-flow using a 60-minute gradient. The samples consisted of three separate proteomes (HeLa, yeast, and *E. coli*), which are mixed in various proportions.

Results

The quality and breadth of the spectral library has a critical impact on the detection of peptides and proteins from DIA data sets. Therefore an exhaustive spectral library is created using multiple off-line high pH reverse phase fractions of each proteome. Quantification by DIA is based on the precursor level in the high resolution full scan. Scan cycles consisting of multiple full- MS1 scans with interspersed DIA scans are employed to improve the precision of MS1 based quantification. For robust and reproducible chromatographic separation, a 60-minute capillary LC- DIAMS method is developed without compromising sensitivity along with high-throughput. The median CV of all the experiments are < 10%. Aroll-up statistics strategy is applied to calculate peptide p-Value, improving quantitation precision to 95% by removing the interferences from the real quantitative signals. As a result, the ratios of the mixed three proteomes accurately reflect the expected values.

Conclusions

An HRMS1–DIA workflow solution has been developed by streamlining each key component to achieve ease-of-operation and reproducible label-free quantitative performance in a high-throughput manner. The initial data demonstrate >7000 proteins were reliably quantified. For a higher proteome coverage, a 120-minute DIAMS method based on the nano-LC separation is developed, where ~20% more protein groups (> 8000) were quantified.

Novel Aspect

With the novel and complete HRMS1-DIA workflow ca. 10,000 protein groups are identified with 1% FDR and ca. 8000 protein groups are quantified with 95% quantitation precision.

MP-194 / ADVANCING MASS SPECTROMETRY-BASED LARGE-COHORT PROTEOMICS FOR PRECISION MEDICINE – AN INTERNATIONAL CANCER MOONSHOT MULTI-SITE STUDY

Yue Xuan⁽¹⁾ - ***Nicholas W. Bateman***⁽²⁾ - ***Yuju Chen***⁽³⁾ - ***Albert Sickmann***⁽⁴⁾ - ***Bernd Wollscheid***⁽⁵⁾ - ***Connie R. Jimenez***⁽⁶⁾ - ***Martin Røssel Larsen***⁽⁷⁾ - ***Hu Zhou***⁽⁸⁾ - ***Siqi Liu***⁽⁹⁾ - ***Zhinan Chen***⁽¹⁰⁾ - ***Thomas Kislinger***⁽¹¹⁾ - ***Ben Crossett***⁽¹²⁾ - ***Sebastien Gallien***⁽¹³⁾ - ***Yue Zhou***⁽¹⁴⁾ - ***Pedro Navarro***⁽¹⁵⁾ - ***Reta Birhanu Kitata***⁽³⁾ - ***Christin Lorenz***⁽⁴⁾ - ***Sandra Goetze***⁽⁵⁾ - ***Sander Piersma***⁽⁶⁾ - ***Davide Chiasserini***⁽⁶⁾ - ***Muhammad Tahir***⁽⁷⁾ - ***Hongwen Zhu***⁽⁸⁾ - ***Guixue Hou***⁽⁹⁾ - ***Xiuxuan Sun***⁽¹⁰⁾ - ***Andrew Macklin***⁽¹¹⁾ - ***Amanda Khoo***⁽¹¹⁾ - ***Benjamin L. Parker***⁽¹⁶⁾ - ***Stuart J. Cordwell***⁽¹⁶⁾ - ***Thomas P. Conrads***⁽²⁾

Thermo Fisher Scientific, CMD, Bremen (1) - Uniformed Services University of the Health Sciences, Gynecologic Cancer Center of Excellence, Murtha Cancer Center, Bethesda (2) - Academia Sinica, Institute of Chemistry, Taipei (3) - Leibniz-Institut für Analytische Wissenschaften – ISAS – e. V, ISAS - Dortmund, Dortmund (4) - Eidgenössische Technische Hochschule Zürich, Institute of Molecular Systems Biology (IMSB), Zurich (5) - VU University Medical Center, Dept. Medical Oncology, Cancer Center Amsterdam, Amsterdam (6) - University of Southern Denmark, Department of Biochemistry and Molecular Biology, Odense (7) - Chinese Academy of Sciences, Shanghai Institute of Materia Medica, Shanghai (8) - BGI-SHENZHEN, BGI-SHENZHEN, Shenzhen (9) - The Fourth Military Medical University, Cell Engineering Research Center / National Key Discipline of Cell Biology, Xi'an (10) - Princess Margaret Cancer Centre, Princess Margaret Cancer Centre, Toronto (11) - The University of Sydney, Sydney Mass Spectrometry,, Sydney (12) - Thermo Fisher Scientific, CMD, Paris (13) - Thermo Fisher Scientific, CMD, Shanghai (14) - Thermo Fisher Scientific, CMD, Bremen (15) - The University of Sydney, School of Life and Environmental Science, Sydney (16)

Keywords: Cancer moonshot multi-site study, Large-cohort study, High Resolution MS1 based Data-independent acquisition (HRMS1-DIA), High-throughput, Reproducibility

Introduction

A high-throughput and streamlined analytical workflow using high resolution MS1-based quantitative data-independent acquisition (HRMS1-DIA) is standardized with well-defined experimental steps, and systematically applied to a set of samples. The study was benchmarked across multiple Cancer Moonshot sites worldwide utilizing same instrument platform, procedures, and software, and demonstrated to be stable in a 24/7 operation mode for 7 consecutive days.

Methods

Experiments were performed on a capillary chromatographic separation with 1h gradient online coupled to Q Exactive HF mass spectrometer. HRMS1-DIA method contains multiple MS1 scans with interspersed DIA MS2 scans. The samples consisted of three separate proteomes mixed in various proportions. HeLa digest served as a quality control (QC) sample. Spectral libraries were created by performing LC-MS/MS analysis of multiple off-line high pH reverse phase fractions.

Results

>250 QC files and >360 mixed proteome sample files are acquired by the participated laboratories in a 24/7 operation mode for 7 consecutive days.

With the 1 hour capillary LC-HRMS1 DIA workflow, > 5,000 protein groups from > 40,000 peptides of the HELA QC sample are consistently identified with 1% FDR across all sites and all 7 days. The peptide intensity correlation of all the 250 QC files is > 0.85, demonstrating an ultrahigh inter-laboratory and inter-day reproducibility.

> 7,000 protein groups from > 50,000 peptides of the mixed proteome sample are reliably quantified across all 11 international labs and 7 days. >80% proteins are overlapped identified and quantified across the laboratories. Additional to the 1% FDR, a rollup strategy is developed and applied to remove the none-reliable peptide quantitation ratios, to successfully enhance the quantitation precision to > 90%. As a result, the label-free quantitation ratios of the mixed three proteomes accurately reflect the expected values at each laboratory across all 7 days.

Conclusions

To elevate discovery proteomics to translational research in the pipeline of precision medicine, large-cohort studies are essential in discovery and verification of protein biomarkers. Yet, reproducible proteome profiling and quantitation in different laboratories remains a challenge. This multi-site study across the globe demonstrates the high-throughput, reproducibility, and scalability of the standardized HRMS1-DIA workflow for large-cohort studies.

Novel Aspect

Multi-site reproducibility, robustness, and method transfer evaluation for HRAM MS-based global proteome profiling using a standardized DIA analytical workflow

MP-195 / GLOBAL ANALYSIS OF PROTEIN COMPLEXES IN TRYPANOSOMA BRUCEI USING CROSS-LINKING MASS SPECTROMETRY

Fernando Almeida - Michele Tinti - Michael Ferguson

The University of Dundee, Biological Chemistry and Drug Discovery, Dundee, United Kingdom

Keywords: protein correlation profiling; Trypanosoma brucei; cross-linking; size exclusion chromatography; cryomilling

Introduction

Trypanosoma brucei is a protozoan parasite and the causative agent of sleeping sickness in humans.

Approximately 9100 protein coding genes have been identified by whole genome sequencing[1], yet 4900 of these lack classifiable homology to known proteins in other organisms. To further ascribe functional classification to these “hypothetical” proteins [2,3], we performed a protein correlation profiling (PCP) based on a cross-linking mass spectrometry approach.

Methods

Cells were cryomilled and the lysate powder was reacted with the cross-linking reagent DTSSP. After, the cross-linked protein samples were submitted to size exclusion chromatography (SEC) followed by their fractionation. Each fraction was submitted to nanoLC-MSMS analyses and their PCP analysis was performed based on label free quantification (LFQ) intensities processed by MaxQuant.

Results

PCP identifies co-fractionating proteins to enable a global characterization of functionally related proteins via “guilt by association”, as many intracellular biological processes are dependent on the stable physical association between two or more proteins. We utilised cryomilling to maintaining vivo protein associations, whilst using chemical crosslinking to stabilises such interactions during further sample processing.

In order to optimise DTSSP mediated crosslinking conditions, α -tubulin was monitored as a reference protein by Western blotting with anti- α tubulin. The established reaction condition of 4 mM DTSSP for 30 mins was sufficient to induce crosslinking. As a next step, the SEC was performed, and the eluted protein complexes were collected in 12 fractions. The PCP results revealed protein complexes that were not seen in the previous study carried out using native conditions[4]. Moreover, we have found hypothetical proteins associated with annotated proteins, that might help to accredit their shared function.

Conclusions

We established a crosslinking approach with DTSSP on cryomilled *Trypanosoma brucei* cell lysates, followed by their proteins separation on SEC and fractionation. Using a small number of fractions, PCP preliminary data have shown the capability to reveal not only protein complexes that were not seen in the previous study, but also the protein associations between annotated and hypothetical proteins.

Novel Aspect

The novel aspect of this approach relies on the combination of cryomilled cell lysates, the covalent protein association using crosslinking, SEC and PCP.

References

1. Berriman, M et al. Science, 309(5733), 416-422 (2005)
2. Salavati, R., and Najafabadi, H. S. Trends Parasitol. 26, 225-229 (2011)
3. Alsford, S. et al. Genome Res. 21, 915-924 (2011)
4. Crozier, T. W. M. et al. Mol. Cell Proteomics, 2254–2267 (2017)

MP-196 / ABUNDANT PROTEIN DEPLETION OF HUMAN PLASMA SAMPLES – SAMPLE PREPARATION APPROACHES FOR QUANTITATIVE COMPARISON STUDIES

Ryan Bomgardner - Sergei Snovida - Katherine Herting - Ramesh Ganapathy - Barbara Kaboord - John C. Rogers

Thermo Fisher Scientific, R&D, Rockford, United States

Keywords: plasma, depletion, TMT, sample preparation

Introduction

Depletion of abundant plasma proteins is required to identify and measure changes in prognostic or

diagnostic proteins [1]. We have developed new top2 and top12 abundant protein depletion resins and evaluated them for depletion specificity, efficiency, and reproducibility. We also combined our abundant protein depletion resins with a streamlined Tandem Mass Tag (TMT) sample preparation procedure.

Methods

Commercially-obtained pooled human plasma samples were used to assess depletion resins by SDS gel, ELSIA, and LC-MS. High pH reversed-phase fractionation was performed on TMT-labeled samples before LC-MS analysis. All samples were analyzed on Thermo Scientific Orbitrap Fusion Tribrid mass spectrometer and processed using Thermo Scientific Proteome Discoverer 2.2 software [2].

Results

After optimization of the antibody ligand conjugation chemistries and resin blending protocols, we achieved >97-99% depletion efficiency of twelve high abundance protein targets in human plasma samples. Abundant protein depletion from plasma enabled detection of 34% more peptides and 25% more proteins compared to undepleted samples. Assessment of depletion capacity of the new resins showed high reproducibility for samples ranging from 10-100ul. These depletion resins were also used in combination with a new sample preparation workflow for TMT-reagent labeling which minimizes sample handling and the number of steps. Using this workflow, we identified and quantified nearly 1100 proteins in normal and lung cancer patient plasma samples.

Conclusions

Abundant protein depletion from plasma enables detection of more peptides and proteins for better detection and quantitation of relevant biomarkers. Combining abundant protein depletion with a streamlined TMT reagent sample preparation workflow and high pH reversed phase fractionation resulted in quantification of nearly 1100 proteins in normal and lung cancer patient plasma samples.

Novel Aspect

Development and comprehensive analysis of new top2 and top12 depletion resins for abundant human plasma proteins for quantitative comparison studies.

References

1. Anderson N. L., Anderson N. G. Mol Cell Proteomics. 11:845-67 (2002).
2. Erickson B., Jedrychowski M., McAlister G., Everley R., Kunz R., Gygi S. Anal. Chem. 87:1241-1249 (2015).

MP-197 / QUANTIFICATION OF HUMAN INSULIN AND ITS ANALOGUES BY IMMUNOPURIFICATION AND LC-MS/MS

Jeroen Kooistra - Josipa Grzetic Martens - Younes Lamkadam - Theo Noij

Charles River Laboratories, Bioanalysis, Den Bosch, Schaijk, Netherlands

Keywords: human insulin, aspart, quantification, LC/MS, immunopurification

Introduction

Diabetes, a disease altering the normal insulin production, affects 422 million of people globally and research into its mechanisms and treatment is a growing research area (1). Current quantifications for insulin and its analogues, such as ELISA, have flaws in terms of performance and throughput (2). Liquid chromatography–mass spectrometry (LC-MS/MS) approaches have recently been developed as alternative strategies for analysis of biomolecules, including insulin.

Methods

Insulin immunocapture followed by LC-MS/MS.

Sample preparation includes capturing of insulin (aspart and human) using immobilized (on magnetic beads) biotinylated antibody from insulin spiked (rat or human) plasma. Insulin is then eluted and subjected to the Glu-C digestion prior LC-MS/MS analysis using Waters UPLC Peptide CSH C18 combined with Sciex Qtrap6500+.

Results: Method validation for quantification of insulin aspart and human insulin

All investigated calibration curves were characterized by a high correlation coefficient ($R^2 \geq 0.99$).

For Insulin aspart a method was developed from 1.75 to 700 ng/mL. Good selectivity of the method was confirmed with only 2.6% relative analyte response in the blank sample. Quality control samples of 1.75, 5.25, 36.8 and 525 ng/mL were analyzed (in fivefold for each level), with high accuracy (97%, 97%, 101% and 95 % respectively), demonstrating good precision of the method (CV=4%, 6%, 11% and 22 % respectively).

The assay for human Insulin quantification is under development and currently viable for quantification of human insulin in BSA/PBS for 0.35 ng/mL to 15 ng/mL. For this assay we have the ultimate goal to go lower to 70pg/mL in human plasma, which is necessary for insulin quantification assays according to the clinical guidelines. By improving the selectivity of method by antibody and internal standard selection, and by LC-MS parameter optimization, this goal is considered to be feasible, on which the method is currently under development.

Conclusions

With significant reduction of matrix interference by immunoaffinity sample preparation and with sensitivity of Sciex Qtrap 6500+, a low level of quantification of 1.75 ng/mL and 0.35 ng/mL for insulin aspart and human insulin, respectively, is achieved. To set the LLOQ to 0.07ng/mL the blank response has to be reduced. To achieve that, additional test are on the way: to use another more selective antibody and to investigate other selective MS transitions.

Novel Aspect: Outlook of the method

This research investigates the possibilities to overcome inabilities of current insulin quantification of human and animal origin and insulin drug analogues by developing a LC-MS assay instead. At Charles River we go towards the implementation of LC-MS protein quantitation to set a sensitive, selective and reproducible method for the quantification of insulin according to the clinical guidelines that is soon to become a fast, reliable and routine analysis.

References

1. World health organization, www.who.int
2. Miller W. G., Thienpont M. L., Van Uytvanghe K., Clark P. M., Lindstedt P., Nilsson G., and Steffes M. W., *Clinical Chemistry*, 55:5, 1011-1018 (2009)
3. Wiśniewski J. R., *Analytical Chemistry*, 88(10), 5438-5443 (2016), doi: 10.1021/acs.analchem.6b00859.

MP-198 / PARALLEL EVALUATION OF MALDI-TOF AND TANDEM LC-ESI MASS SPECTROMETRY FOR RAPID QUANTIFICATION OF ABRIN TOXIN IN COMPLEX MATRICES

Sandrine Livet⁽¹⁾ - **Sylvia Worbs**⁽²⁾ - **Eva-Maria Hansbauer**⁽¹⁾ - **Hervé Volland**⁽¹⁾ - **Stéphanie Simon**⁽¹⁾ - **Christophe Junot**⁽¹⁾ - **François Fenaille**⁽¹⁾ - **Brigitte G. Dorner**⁽²⁾ - **François Becher**⁽¹⁾

CEA Saclay, Service de Pharmacologie et Immunoanalyse (SPI) - Laboratoire d'Etude du Métabolisme des Médicaments, Gif-sur-Yvette, France (1) - Robert Koch Institute, Biological Toxins - Centre for Biological Threats and Special Pathogens, Berlin, Germany (2)

Keywords: toxin, MALDI, PRM, proteomics, bottom-up

Introduction

Biological toxins are relevant in the security, health and food sectors. Development of new analytical methods for rapid and unambiguous identification of low toxin concentrations in various complex matrices

is essential. We reported previously a method based on LC-ESI-MS/MS for abrin quantification in food and environmental matrices [1]. Here, immuno-MALDI-TOFMS was evaluated for rapid and chromatography-free quantification of abrin.

Methods

Magnetic beads coated with abrin-specific antibodies were mixed with 500 µL of food samples spiked with the toxin. After harsh washing conditions with PBS-tween, analytes were eluted from beads by incubation with 0.1% TFA. Abrin was heated for 15 min at 95°C in RapiGest SF 0.05% for denaturation and digested with trypsin in a bath-type sonicator. Peptides were co-crystallized with CHCA matrix on a disposable plate for MALDI-TOF analysis using a Bruker UltrafleXtreme.

Results

Up to 29 peptides of abrin isoforms were detected by MALDI-TOF MS in a pure standard and the 4 most intense, different to the best ESI peptides, were selected. Isotope-labeled peptides were added for quantification. Immuno-capture conditions were improved for MALDI-TOF MS, including additional washes of antibody-coated beads. On-beads digestion resulted in high interfering signals from antibody peptides, and abrin was therefore eluted before digestion. Trypsin amount was reduced 10 times and digestion time increased to 1 hour in the bath-type sonicator. Using the new protocol, interfering peaks and ionization suppression in apple juice and milk samples were decreased significantly, resulting in spectra with predominant abrin peptide signals at 250 ng/ml. At this concentration, about twenty peptides of abrin isoforms were detected in the two food matrices. On-going parallel evaluation with LC-ESI-MS/MS shows sensitivity in the same range with limit of quantification in milk at 10 and 40ng/ml for LC-ESI-MS/MS and MALDI-TOF, respectively.

Conclusions

Bottom-up proteomics and immuno-MALDI-TOF MS allowed quantification of abrin and its isoforms in complex matrices, without time consuming peptide chromatography or extraction resulting in simple and fast (< 5 min) MS analysis. The method was found specific, reliable, and sensitive in milk. Further optimizations are on-going to reach equal or better limit of quantification than LC-ESI-MS/MS and to reduce the total analysis time by automatization of the workflow.

Novel Aspect

Evaluation of MALDI-TOF and LC-ESI-MS/MS for rapid and sensitive quantification of abrin in complex matrices.

References

1. Hansbauer EM., Worbs S., Volland H., Simon S, Junot C., Fenaille F., Dorner B. G., and Becher F., Analytical chemistry, 89, 11719-11727 (2017).

MP-199 / PROTEOMIC ANALYSIS OF UNCONVENTIONAL PROTEIN SECRETION

Corina Mayrhofer⁽¹⁾ - ***Katharina Nöbauer***⁽²⁾ - ***Gottfried Brem***⁽³⁾

University of Veterinary Medicine Vienna, Institute of Animal Breeding and Genetics / Institute of Biotechnology in Animal Production, IFA Tulln, Vienna, Austria (1) - University of Veterinary Medicine Vienna, VetCore Facility for Research, Vienna, Austria (2) - University of Veterinary Medicine Vienna, Institute of Animal Breeding and Genetics, Vienna, Austria (3)

Keywords: autophagy-related protein 7, unconventional protein secretion, epithelial ovarian cancer, secretome

Introduction

The cancer secretome comprises numerous proteins that promote tumorigenic processes [1]. A substantial

number of proteins lack signal peptides and are released via unconventional mechanisms[2]. Herein we analyzed the contribution of autophagy-related (ATG) based unconventional protein secretion to the microenvironment of epithelial ovarian cancer cells.

Methods

Epithelial ovarian cancer cells (OV-90, ATCC) were transfected with siRNA targeting human ATG7 or control siRNA. A label free quantitative approach was applied to analyse secreted proteins of ATG7 deficient cells and ATG7 intact (control) cells. Conditioned media was collected over 24 hours. Samples were digested using FASP [3] and analysed by nanoLC-MS/MS (Orbitrap MS). Raw data processing was performed with MaxQuant.

Results

ATG7 siRNA transfection efficiently reduced ATG7 protein levels. The secretome was highly enriched with extracellular localized proteins. Approximately one-third of these proteins contained a signal peptide as predicted by Secretome P. Multivariate analysis of the data further demonstrated changes of the secretory protein profile upon siRNA mediated ATG7 knockdown. A subset of proteins emerged that were not quantified in the conditioned media of ATG7 deficient cells, including proteins involved in cell-cell adhesion.

Conclusions

Besides the well-known function of ATG7 in the elongation reaction and the formation of the autophagosome [4], the present proteomic approach shows the contribution of ATG7 protein to the extracellular release of a subset of proteins in epithelial ovarian cancer cells.

Novel Aspect

Novel insights into the complexity of the epithelial ovarian cancer cell microenvironment.

References

1. Vaughan, S., Coward, J. I., Bast, R. C., Jr., Berchuck, A., et al. *Nat Rev Cancer*, 11, 719-725 (2011).
2. Rabouille, C., Malhotra, V., Nickel, W. *J Cell Sci*, 125, 5251-5255 (2012).
3. Wisniewski, J. R., Zougman, A., Nagaraj, N., Mann, M. *Nat Methods*, 6, 359-362 (2009).
4. Green, Douglas R., Levine, B. *Cell*, 157, 65-75 (2014).

MP-200 / PROTEOME OF AVIAN PATHOGENIC E. COLI DURING INFECTION: METHODS AND RESULTS

Arkadiusz Nawrocki ⁽¹⁾ - **Andreas Eske Johansen** ⁽²⁾ - **Jakob Møller-Jensen** ⁽¹⁾ - **John Elmerdahl Olsen** ⁽²⁾ - **Martin Røssel Larsen** ⁽³⁾

University of Southern Denmark, Biochemistry and Molecular Biology, Odense, Denmark (1) - University of Copenhagen, Department of Veterinary and Animal Sciences, Frederiksberg, Denmark (2) - University of Southern Denmark, Biochemistry and Molecular Biology, Odense, Cyprus (3)

Keywords: Proteomics; Mass spectrometry; Bacterial infections;

Introduction

Avian pathogenic E. coli (APEC) gives rise to colibacillosis which result in economic losses in the production system and reduced animal welfare. At the same time, an extensive use of antibiotics to combat infections promotes selection of resistant bacteria. We investigated the proteome of an E. coli ST117 strain during infection of the mucus layer of oviduct to suggest potential genes (essential and redundant) as targets in development of antibiotics or helper drugs.

Methods

A combination of reducing agent and a number of enzymes was used to purified and enriched bacteria

from mucus in quantities enough for LC-MS based quantitative proteome analysis. Peptides from infecting bacteria and the host (in the same sample) were labelled with TMT isobaric tags and mixed with peptides derived from proteins of ST117, grown in culture media (log and stationary phase), labelled with complementary TMT tags.

Results

By combining reducing agent and selected enzymes we successfully purified and enriched bacteria from mucus of infected hens. We identified over 2000 bacterial proteins and were able to quantify approx. 1500 of them across all the samples. This relatively high number of positive hits was greatly facilitated by spiking TMT labelled peptides derived from proteins extracted from ST117 grown in culture media into these obtained from infection experiments. The abundant spiked in TMT labelled peptides served as factor to isolate E. coli derived peptides for fragmentation in the LC-MSMS analysis, as otherwise peptides from co-isolated proteins of the host highly suppressed the signal of peptides originating from isolated bacterial proteins.

Conclusions

The approach presented here is suitable for detection and quantification of relatively high number of proteins of bacteria isolated from infected oviduct of hens, where proteins of the latter highly dominate the sample. Apparently, this approach could be adopted to investigate pathogen/host interactions, or proteomes of mixed population of different organisms.

Novel Aspect

Improved method for isolation of mucus infecting bacteria. Combination of TMT labelled peptide populations derived from pure cell culture and isolates from infected hosts.

MP-201 / COMPARATIVE STUDY OF THREE MS-BASED APPROACHES FOR QUANTITATIVE PROTEOMICS ANALYSIS

Domenico Taverna - Marco Gaspari - Marco Catracchia

Proteomics-UMG lab, Department of Experimental and Clinical Medicine, Magna Graecia University of Catanzaro – ITALY

Keywords: quantitative proteomics, TMT, LFQ, dimethyl labeling

Introduction

Several approaches are available for performing MS-based quantitative proteomics. These include label-free and label-based methods [1]. However, which quantification strategy is best suited for identification of biologically relevant changes at the protein level is still not clear [2]. In this study, the analytical performance of three different strategies for relative protein quantification using shotgun proteomics have been compared: LFQ, DML and TMT.

Methods

Two protein mixtures were analyzed, both containing the same amount of human proteins, but different protein amounts from E. Coli (1.5-fold difference) and soybean (3-fold). After protein denaturation and tryptic digestion (5 technical replicates), peptides were divided in three aliquots: one was used for label-free analysis; the second was labeled by TMT, while the third was derivatized by dimethyl labeling. All samples were analyzed by LC-MS/MS on a Q-Exactive.

Results

Data were processed for statistical analysis (MaxQuant/Perseus) aiming to find out the best conditions for

acceptable diagnostic accuracy. For 5 replicated samples per protein mix, a total of 2495 (LFQ), 1385 (DML), and 1803 (TMT) proteins with minimum 2 peptides (FDR 0.01) were quantified. Parameters such as specificity, precision and sensitivity were evaluated for each of the three quantification strategies. Differentially expressed proteins between protein mix A and B for LFQ and TMT datasets were identified by using a corrected Student t-test (permutation-based FDR < 0.05, S0=0.2). Instead, DML data were processed using the following criteria: one sample t-test, S0=0.2, p-value <0.05. The diagnostic accuracy has been evaluated for all the three datasets: LFQ data returned 95% specificity, 71% precision, and 87% sensitivity; samples labeled with TMT obtained 100% specificity, 96% precision and 51% sensitivity; DML labeling instead gave back the following values: 96% specificity, 64% precision, 79% sensitivity.

Conclusions

Since LFQ and DML provided comparable diagnostic accuracy, in our analytical pipeline LFQ was preferred because it provided deeper proteome coverage and faster sample preparation. TMT provided excellent precision (low number of false positives), though sensitivity was poor, especially at low fold changes (1.5-fold).

Novel Aspect

Comparison of the diagnostic accuracy of the three main categories in quantitative proteomics: label-free, isobaric tagging and label-based with full MS quantification.

References

1. Li H., Han J., Pan J., Liu T., Parker C. E., Borchers C. H., J. Mass Spectrom., 52, 319–341 (2017).
2. Hogrebe A., von Stechow L., Bekker-Jensen D. B., Weinert, B. T., Kelstrup C. D., Olsen J. V., Nat. Commun., 9, 1045 (2018).

MP-202 / PROTEOMIC AND METABOLOMIC ALTERATIONS IN THE RAT BRAIN AFTER ISCHEMIC STROKE

Jaakko Teppo⁽¹⁾ - *Anu Vaikkinen*⁽¹⁾ - *Kert Mätlik*⁽²⁾ - *Tapio Kotiaho*⁽¹⁾ - *Mikko Airavaara*⁽²⁾ – *Risto Kostainen*⁽¹⁾

University of Helsinki, Faculty of Pharmacy, Helsinki, Finland (1) - University of Helsinki, Institute of Biotechnology, Helsinki, Finland (2)

Keywords: ischemic stroke, proteomics, metabolomics

Introduction

Stroke is one of the major causes of death worldwide, and in non-lethal cases often leads to permanent symptoms[1]. Understanding the biological phenomena occurring in the brain after stroke on a molecular level is crucial for the development of new treatments to alleviate these symptoms. The aim of our study was to characterize the effects of stroke on the proteome and metabolome of rat brain at a timepoint in which the brain's response to stroke has been active for days.

Methods

Rats were subjected to 60 min cerebral ischemia or a sham operation as a control. The peri-infarct areas were dissected 4 days later, and the samples were homogenized and divided into two aliquots. Sample preparation was done separately for proteomics and metabolomics, and the untargeted analyses were performed with Thermo EASY-nLC 1000/Q Exactive MS and Dionex UltiMate 3000 UHPLC/Orbitrap Fusion MS, respectively. Data was analyzed with MaxQuant, MZMine2, and R.

Results

Hierarchical clustering and PCA of the data indicated major differences between the sample groups, and in

the case of PCA, higher between-sample variability in the stroke group than in the control group in both proteomics and metabolomics. Across the samples, in total 2403 proteins were quantified, of which 277 were upregulated and 167 downregulated in the stroke samples relative to the control samples ($q < 0.01$). In line with the literature, inflammatory response-related proteins were upregulated in the brain after stroke[2]. Further processes identified by the bioinformatic analysis include *e.g.* increased translation and decreased energy metabolism in the stroke samples. In metabolomics, 470 features were significantly different ($q < 0.01$, fold change > 2) between the sample groups. Based on MS/MS identifications, the abundance of *e.g.*, sugars and sugar phosphates, unsaturated fatty acids and a wide variety of other lipid species increases in the brain after stroke, whereas the abundance of nucleosides (adenosine, guanosine) decreases.

Conclusions

Ischemia causes extensive damage to the brain, and with the consequent reperfusion, induces various biological processes. Our work identified major alterations in the brain after stroke on both proteomic and metabolomic level, including proteins and metabolites related to inflammation, protein synthesis, and energy metabolism.

Novel Aspect

This study integrates proteomics and metabolomics data from the peri-infarct brain area of rats subjected to cerebral ischemia and sampled several days post-stroke.

References

1. Hankey G. J., *The Lancet*, 389, 641-645, 2017
2. Anrather J., Iadecola C., *Neurotherapeutics*, 13, 661-670, 2016

MP-203 / DETERMINATION OF PEPTIDE PURITY BY UV-MS AND QUANTITATIVE NMR

Frank A. Torma - Milena Quaglia - John Warren

LGC LTD, Chemical Measurement and Calibration, Teddinton, United Kingdom

Keywords: peptide, quantification, SI-traceable, UV-MS, standardization

Introduction

Accurate quantification of peptides and proteins rely on the availability of certified reference materials. Here we investigated the use of ultraviolet (UV) and mass spectrometric (MS) detection for the simultaneous determination of the quantity and identity of peptidic impurities present in peptide standard solutions. The data combined with amino acid analysis can provide accurate determination of the SI-traceable concentration of peptide calibrants.

Methods

Peptide solutions were SI-traceably quantified by amino acid analysis using exact matching isotope dilution mass spectrometry (EM-IDMS). The relative amount and the identity of the impurities was determined by UV-MS. Relative response factors (RRF) were verified by varying the injection volumes of the stock solution. The amino acid results were corrected for peptidic impurities. The results were confirmed by quantitative NMR (qNMR).

Results

The SI-traceable concentration of the peptide solutions can be determined by peptidic impurity corrected amino acid analysis (PICCA). The peptide and its impurities must absorb UV light and the impurities must be baseline-separated. The peptidic impurities were sequenced by MS/MS. The main impurities are deletion

products and related sequences with post translational modifications on unstable amino acid residues. The determined purity of the synthetic peptides used in this study is indirectly proportional with the length of the peptide sequence. The validation of the method for the determination of the purity of potential peptide calibrators will take place by applying it to the determination of the purity of oxytocin as part of an international comparison study. Quantitative NMR results complement the findings of the UV-MS results.

Conclusions

The relative (detector) response ratios for the impurities vs the main compound can be determined without the use of solid impurity standards, and the determination of the purity of a peptide solution by the evaluation of the UV chromatogram is possible. The agreement between the purity results determined by UV and MS using the orbitrap mass analyzer suggests linear behavior for the orbitrap.

Novel Aspect

Determination of purity of synthetic peptides is achieved with the use of universal (UV) and mass spectrometric detection and quantitative NMR.

MP-204 / UNIVERSAL S-TRAP SAMPLE PROCESSING: STANDARDIZED REPRODUCIBLE SAMPLE WORKUP FROM < UG TO > MG SCALES FOR ALL SAMPLE TYPES

John Wilson - Darryl J. Pappin

Cold Spring Harbor Laboratory, Cold Spring Harbor, NY, USA

Keywords: sample preparation, digestion, sample cleanup, bottom-up proteomics, FFPE

Introduction

Sample preparation has long hampered bottom-up proteomics. Proteins' vastly different solubilities and myriad extraction protocols yield very different data depending on what proteins are [not] extracted. S-Trap sample processing[1] solves this problem by integrating: 1) complete 5% SDS protein solubilization; 2) total protein denaturation (pH < 1 and > 70% organic); 3) simultaneous sample concentration/cleanup; and 4) rapid reactor-type proteolytic digestion.

Methods

S-Traps were obtained from ProtiFi (www.protifi.com, Huntington NY). S-Trap micros were used from sub-ug to 100 ug max, S-Trap minis from 100 – 300 ug, S-Trap midis > 300ug, and 96 well plates for high throughput applications. Briefly, all samples were extracted with 5% SDS, fully denatured, captured on the S-Trap where they were cleaned of all contaminants including Laemmli sample buffer, PEG, detergents, urea, reduction and alkylation reagents, etc. and digested.

Results

S-Trap sample processing allowed unbiased protein analysis including poorly soluble proteins(membrane proteins, chromatin and muscle, etc.) and was successfully applied without alteration to all samples including tissues, FFPE blocks, cell lines, biological fluids, etc. Using only standard lab equipment, S-Traps enabled reproducible sample processing with < 10% CVs from protein to peptides in around 1 – 2 hrs from sub-microgram to multi-milligram scales.

Conclusions

S-Trap processing is highly reproducible (CVs typically < 10%) and independent of the properties of the proteins under study. Without change, its harmonized protocol is applicable to all kinds of samples from liquid biopsies (serum and dried blood spots [DBS]) to FFPE blocks to tissues or cell lines. S-Trap processing increases protein recovery, affording analysis of very small sample quantities, and can enable simultaneous proteomics and metabolomics.

Novel Aspect

Use one protocol for all samples from <ug to >mg scales greatly improves reproducibility and ease of research. Throughput is greatly improved by 1 hr digests. The ability to remove contaminants is especially helpful to core labs.

References

1. Zougman A, Selby PJ, Banks RE. *Proteomics*. 2014 May;14(9):1006-0.

MPS-S08 - Chiral discrimination

Chairs: Anne Zehnacker-Rentien, Maurizio Speranza

MP-205 / SELF-ASSEMBLED BINUCLEAR CU(II)–HISTIDINE COMPLEX FOR ABSOLUTE CONFIGURATION AND ENANTIOMERIC EXCESS DETERMINATION OF NAPROXEN BY TANDEM MASS SPECTROMETRY

Xiangying Yu⁽¹⁾ - *Man-Chu Chau*⁽²⁾ - *Wai Kit Tang*⁽²⁾ - *Chi-Kit Siu*⁽²⁾ - *Zhong-Ping Yao*⁽³⁾

School of Chemical Engineering and Energy Technology, Dongguan University of Technology, Dongguan (1) - Department of Chemistry, City University of Hong Kong, Kowloon Tong (2) - Department of Applied Biology and Chemical Technology, The Hong Kong Polytechnic University, Hung Hom (3)

Keywords: Binuclear Cu(II)-L-Histidine Complex, S/R-Naproxen, Chiral Differentiation, Ion-Mobility Spectrometry, DFT Molecular Dynamics Modeling

Introduction

Naproxen is one of the most consumed nonsteroidal anti-inflammatory drugs and marketed as S-naproxen since its R-enantiomer is hepatotoxic[1]. Developments of efficient and reliable methods for chiral differentiation and quantification of naproxen are thus essential [2].

Methods

Recognition and quantification of S/R-naproxen enantiomers with commonly used chiral selectors are studied using tandem mass spectrometry (MS/MS). Structures of the most promising diastereomeric precursors are revealed using collision-induced dissociation (CID), ion-mobility spectrometry (IMS), and density functional theory (DFT) and molecular dynamics(MD).

Results

Among all studied diastereomeric complexes formed between S/R-naproxen and various chiral selectors [3], including $\alpha/\beta/\gamma$ -cyclodextrins, amino acids and their derivatives, glucose, tartaric acid, and vancomycin, a novel binuclear metal bound $[(M^{II}\text{-His})_2(\text{S/R-naproxen})-3\text{H}]^+$ complex (M = Cu, Ni, or Co with Cu being the best) can allow effective absolute configuration identification and enantiomeric excess determination for naproxen using MS/MS. The key candidate structure of $[(\text{Cu}^{II}\text{-His})_2(\text{S/R-naproxen})-3\text{H}]^+$ revealed from CID, IMS and DFT-MD, has an unusual self-assembled, compact geometry with the two Cu(II) ions bridged closely together by the carboxylate groups of the two histidines. The difference in dissociation efficiency of the two diastereomers is attributed to the σ - π interaction between the NH_2 group of one histidine and the naphthyl ring of naproxen in the complex.

Conclusions

The present report observes and characterizes the diastereomeric complexes of $(\text{Cu}^{II}\text{His})_2$ with aromatic acid, which could contribute to the recognition of chiral aromatic acids, design of catalysts based on

binuclear copper-bound complex, as well as the better understanding of metal-ion complexation by His or His-containing ligands.

Novel Aspect

A novel self-assembled, compact $[(\text{Cu}^{\text{II}}\text{L-His})_2(\text{S/R-naproxen})-3\text{H}]^+$ complexes are used for efficient chiral differentiation of naproxen.

References

1. C. Caballo, M. D. Sicilia, S. Rubio, *Anal. Bioanal. Chem.* 407, 4721–4731 (2015).
2. D. Camacho-Munoz, B. Kasprzyk-Hordern, *Anal. Bioanal. Chem.* 407, 9085–9104 (2015).
3. X. Yu, Z. Yao, *Anal. Chim. Acta* 968, 1–20 (2017).

MP-206 / SPECTROSCOPIC DISCRIMINATION OF DIASTEREOMERIC COMPLEXES INVOLVING AN AXIALLY CHIRAL RECEPTOR

Antonello Filippi⁽¹⁾ - **Caterina Fraschetti**⁽¹⁾ - **Laura Guarcini**⁽¹⁾ - **Costantino Zazza**⁽²⁾ - **Tadashi Ema**⁽³⁾ - **Maurizio Speranza**⁽¹⁾

Sapienza-University of Rome, Chemistry and technology of drugs department, Rome, Italy (1) - Sapienza-University of Rome, Chemistry Department, Rome, Italy (2) - Okayama University, Graduate School of Natural Science and Technology, Okayama, Japan (3)

Keywords: IRMPD, chiral discrimination, axial chirality, noncovalent interactions

Introduction

Different stability and reactivity of diastereomeric aggregates reflect the number and intensity of the interactions between the chiral species involved. Therefore, investigating diastereomeric aggregates in the gas phase helps to understand the mechanism of chiral recognition through detailed knowledge of the structure of the chiral receptor, as well as of the size- and shape-specific non-covalent interactions in its aggregate with a chiral molecule [1-3].

Methods

An ESI-MS-IRMPD analysis has been employed to assess the structural features of diastereomeric complexes [4] between an axially chiral receptor (M) and two amino acids, phenylalanine (Phe) and dopa (DOPA), differing only for the aromatic substituents. The results have been interpreted through a computational analysis at the B3LYP/6-31G** level, to assess both the coordination site and the preferential conformation adopted by the complexed guest.

Results

IRMPD spectra of the $[\text{M}\cdot\text{H}\cdot\text{Phe}]^+$ diastereomeric complexes did not exhibit any appreciable difference, in contrast to the significant discrimination exerted by M on the enantiomers of dopa. The comparison between experimental bands and theoretical results indicate that the guest is located out of the cavity, in close proximity of the binaphthyl moiety of M. Furthermore, the computational analysis performed on several rotamers pointed out that the enantiodiscrimination is mostly due to the short range interaction between the hydroxyl groups and the core of the axial chirality. It is worth noting that the structures spectroscopically observed are not the most stable ones.

Conclusions

The IRMPD of the ESI-formed $[\text{M}\cdot\text{H}\cdot\text{DOPA}]^+$ diastereoisomers were consistent with an out of cavity structure, not corresponding to the lower energy coordination site. The observed discrimination on the dopa enantiomers is due to the hydroxyl groups, absent in the structure of phenylalanine, that are oriented

towards the chiral core of M. Finally, the different spectroscopic signatures indicate the coexistence of several conformers for the diastereomeric aggregates.

Novel Aspect

Subtle structural differences between non-covalent diastereomeric complexes have been pointed out through IRMPD technique.

References

1. Schalley C. A., *Mass Spectrom. Rev.*, 20, 253 (2001).
2. Speranza M., Gasparrini F., Botta B., Villani C., Subissati D., Fraschetti C., Subrizi F., *Chirality*, 21, 69 (2009).
3. Zehnacker A., *Chiral recognition in the gas phase*, CRC Press, Boca Raton, FL (2010).
4. Filippi A., Fraschetti C., Piccirillo S., Rondino F., Botta B., D'Acquarica I., Calcaterra A., Speranza M., *Chem. Eur. J.*, 18, 8320 (2012).

MP-207 / ACID - BASED RECOGNITION OF DIASTEREOISOMERIC PEPTIDES USING ION MOBILITY MASS SPECTROMETRY

Magdalena Zimnicka

Polish Academy of Sciences, Institute of Organic Chemistry, Warsaw, Poland

Keywords: ion mobility mass spectrometry, peptides, acids, diastereoisomers

Introduction

The chiral molecules frequently remain undistinguishable using ion mobility mass spectrometry (IM-MS), due to insufficient difference of their collision cross sections at an available mobility resolution of ion mobility drift tubes. The influence of the complexation with organic acids on the ion mobility separation of selected peptide diastereoisomers was evaluated.

Methods

IM-MS measurements were performed on a commercially available (Waters) quadrupole traveling-wave ion mobility time-of-flight spectrometer (Synapt G2-S HDMS). The initial structures of noncovalent complexes obtained with Metropolis Monte Carlo simulations were further optimized at the DFT level of theory. The theoretical cross sections of the optimized associates were calculated using the trajectory method (TM) in the open source software program MOBCAL.

Results

Peptide diastereoisomers containing arginine residue form stable associates with various organic acids in the gas phase. These associates are also stable in conditions of increased pressure of the ion mobility drift cell. Overall the increase of the diastereoisomeric peptide separation after complexation with selected acids was observed. The structures, selectivity, specificity, and the overall effect of complexation with acids on ion-mobility separation of the analyzed diastereoisomers are discussed. Additional relative gas-phase stability experiments were performed to define the association mode between examined peptides and selected organic acids.

Conclusions

The increase of the diastereoisomeric peptide separation after complexation with selected acids was observed. The gas-phase stability studies revealed that the presence of the arginine in a peptide sequence ensures the formation of stable complex in the gas phase. The presence of the available backbone hydrogens in peptide sequence additionally increases the stability of the complexes in the case of diacids.

Novel Aspect

Complexation with the selective acids induces improvement of chiral separation of peptides.

MPS-S09 - Ambient MS and new ionization Methods

Chairs: Jentaie Shiea, Pietro Traldi

MP-208 / A VERSATILE INTEGRATED AMBIENT IONIZATION SOURCE (IAMIS) PLATFORM AND ITS APPLICATIONS

Yu Bai - Wanpeng Ai - Shuting Xu - Huwei Liu

Peking University, College of Chemistry, Beijing

Keywords: ambient ionization, hybrid ionization mode, complex sample, sensitive detection

Introduction

Ambient mass spectrometry (AMS) has become more and more attractive because of its simple sample pretreatment process, the ability to operate in the open air, and capability to probe the surface of samples of any size and shape. The development of AMS opens a new insight for realizing high-throughput, nondestructive, in situ sample analysis from complicate sample[1-3].

The pursuit of high-throughput sample analysis from complex matrix makes the development of multiple ionization techniques with the complementary specialties on demand. Compared with single ambient ionization, combined or hybrid ambient ion source have been attempted to achieve complementary and high-throughput analysis, the reduction of sample complexity, and allowing the simultaneous detection of a broader range of analytes.

Methods

A versatile integrated ambient ionization source (iAmIS) platform was proposed in this work, based on the ideas of integrating multiple functions, enhancing the efficiency of current ionization techniques, extending the application and decreasing the cost of the instrument. The design of the iAmIS platform combines 5 different ionization mode. All individual and combined ionization mode can be easily realized by modulating parameters.

Results

AMS and their combination with other fast sample pretreatment approaches (for example, SPE and SPME etc) and the coupling with different instrument have been developed in our lab [2-4]. By using this iAmIS platform, coconut essential oil and perfume sample has been successfully analyzed. On the other hand, together with the immunal recognition, DESI mode and signal amplification strategy, ultrasensitive detection of proteins in attomolar or zeptomolar were realized.

Conclusions

A versatile integrated ambient ionization source (iAmIS) platform was developed. Complex sample can be analyzed by using this platform, while the sensitive detection of protein at attomolar level can be realized using DESI-MS mode.

Novel Aspect

A versatile integrated ambient ionization source (iAmIS) platform was developed for diverse applications.

References

1. Cooks R. G., Ouyang Z., Takats Z., Wiseman J. M., *Science*, 311:1566-1570 (2006).
2. Zhang J. L., Zhou Z. G., Yang J. W., Zhang W., Bai Y., Liu H. W., *Anal. Chem.* 84: 1496-1503 (2012).
3. Zhang Y. D., Li X. J., Nie H. G., Yang L., Li Z., Bai Y., Niu, L., Song D. Q., Liu H. W., *Anal Chem*, 87: 6505–6509 (2015).

This work was financially supported by the National Natural Science Foundation of China (21575007, 21527809 and 21322505)

MP-209 / THE ROLE OF NEW MATERIALS AS SAMPLING SUBSTRATES FOR THE DESI-HRMS DETERMINATION OF NEW PSYCHOACTIVE SUBSTANCES IN SALIVA

Federica Bianchi⁽¹⁾ - **Silvia Agazzi**⁽¹⁾ - **Nicolò Riboni**⁽¹⁾ - **Minna Hakkarainen**⁽²⁾ - **Leopold Ilag**⁽³⁾ - **Rossana Cecchi**⁽⁴⁾ - **Luca Anzillotti**⁽⁴⁾ - **Maria Careri**⁽¹⁾

University of Parma, Department of Chemistry, Life Sciences and Environmental Sustainability/CIPACK, Parma (1) - KTH Royal Institute of Technology, Department of Fibre and Polymer Technology/Division of Polymer Technology, Stockholm (2) - Stockholm University, Department of Environmental Science and Analytical Chemistry, Stockholm (3) - University of Parma, Department of Medicine and Surgery, Parma (4)

Keywords: Desorption electrospray ionization-high resolution mass spectrometry, sampling substrates, new psychoactive substances, screening method

Introduction

New psychoactive substances (NPS) are a large group of drugs of abuse: being considered as a major threat to public health, the development of rapid methods for their determination is in great demand. DESI-MS enables the in situ analysis of specimens in their native ambient state with minimal sample pretreatment [1,2]. In this context, the role of new materials as sampling substrates is of paramount importance to increase the instrumental response of the analytes.

Methods

Different sampling substrates, *i.e.* commercial PTFE, poly(lactic acid) (PLA)- and silica-based films were synthesized and tested for the DESI-HRMS analysis of NPS (synthetic cannabinoids, synthetic cathinones, ketamine) in saliva. Oral fluid samples from volunteers outside clubs and during roadside patrol control on drivers were collected. A Micro Extraction by Packed Sorbent (MEPS)-DESI-HRMS method was optimized and validated for NPS determination.

Results

The first phase of the study was devoted to the optimization of the DESI parameters affecting the ionization yield and the ion recovery in terms of: i) spray composition; ii) spray flow; iii) acquisition mode; iv) geometrical parameters. With the aim of achieving high sensitivity, different sampling substrates, *i.e.* unfunctionalized PLA, carbon/PLA-based films and a silica-based coating were synthesized and tested. Both unfunctionalized PLA and the silica-based coating proved to be the best choice for NPS detection. A full factorial design followed by the multi-criteria method of desirability functions was applied for the optimization of the MEPS extraction in terms of both loading and eluting cycles. The method was validated obtaining detection limits in the µg/kg level, good linearity, precision and accuracy. Finally, the method was applied for NPS determination in saliva samples: achieving results in agreement with those obtained by a MEPS-GC-MS method used for confirmatory purposes.

Conclusions

The synthesized support materials showed enhanced performances in improving the DESI-HRMS responses

of the investigated compounds. Both the MEPS-DESI-HRMS and MEPS-GC-MS methods represent suitable tools for NPS analysis in oral fluids, combining screening and confirmatory purposes. In particular, the DESI-HRMS method allowed high throughput analyses characterized by both high sensitivity and selectivity, very short analysis times and low costs.

Novel Aspect

Pivotal role of new supporting materials in enhancing the MS responses of the investigated compounds, thus allowing their rapid detection at trace levels in complex matrices

References

1. Takàts Z., Wiseman J. M., Gologan B., Cooks R. G., *Science*, 306, 471 (2004)
2. Bianchi F., Gregori A., Braun G., Crescenzi C., Careri M., *Anal. Bioanal. Chem.*, 407,931 (2015)

MP-210 / MATRIX EFFECTS IN LOW TEMPERATURE PLASMA IONIZATION (LTPI) MASS SPECTROMETRY

Andreas Kiontke⁽¹⁾ - **Birte Scholz**⁽²⁾ - **Walther Honscha**⁽²⁾ - **Claudia Birkemeyer**⁽¹⁾

University Leipzig, Institute of Analytical Chemistry, Leipzig, Germany (1) - University Leipzig, Institute of Pharmacology, Pharmacy and Toxicology, Leipzig, Germany (2)

Keywords: low-temperature plasma ionization, ambient mass spectrometry, matrix effects

Introduction

In recent years, various ionization techniques for mass spectrometry (MS) have been published that allow fast and simple analysis of samples under ambient conditions. Plasma-based methods in particular aroused the interest of many researchers because they offer great potential for use in mobile analytical devices. However, systematic investigations of signal response with this technique, such as the influence of the matrix on signal intensity, are still rare.

Methods

In a first approach, six anilines each in nine different solvents were spotted on a paper target and analyzed by low temperature plasma ionization (LTPI) MS under optimized conditions [1]. In a second approach, 1 µL of a 1 mM chlorpyrifos solution, a common pesticide, was analyzed from peels of different fruits, vegetables and mushrooms. The samples were analyzed using a plasma source coupled to a Bruker 3000+ ion trap.

Results

Our results show that the impact of compound characteristics strongly outperforms the influence of simple sample matrices such as different organic solvents and water, with the trend that volatile solvents tend to decrease the signal responsiveness of the analytes. However, several specific solvent-analyte interactions occurred, which have to be considered in targeted applications of this method.

When choosing a single compound to assess matrix effects, the pesticide chlorpyrifos was detected on all tested different fruits, vegetables and mushrooms. However, we found a strong influence of the matrix on the signal response. For example, citrus fruits showed a particularly strong signal suppression.

Conclusions

For the influence of the solvents, with the exception of water, most solvents seemed to be almost equally well suited for mass spectrometry with LTPI. However, a better signal intensity of the analyte with higher boiling solvents was achieved indicating that a high boiling matrix may be optimal with this ionization technique. In contrast, when investigating the influence of biological matrices such as peels of different fruits, vegetables and fungi, the response of chlorpyrifos was particularly decreased when analysis was

done from the waxy surfaces of citrus fruits. To find the reasons for this behavior is an ongoing project in our lab.

Novel Aspect

We present new insights on matrix effects in low temperature plasma ionization for mass spectrometry.

References

1. Kiontke A., Holzer F., Belder D., Birkemeyer C., *Analytical and Bioanalytical Chemistry*, (2018). DOI: 10.1007/s00216-018-1033-7

MP-211 / APPLICATION OF LDTD/TQMS-HRMS METHODS FOR QUANTIFICATION OF PSYCHOACTIVE COMPOUNDS IN TISSUE MICROSAMPLES

Adam Borik - Andrea Vojs Staňová - Katerina Grabicova - Tomas Randak - Roman Grabic

University of South Bohemia in Ceske Budejovice, Faculty of Fisheries and Protection of Waters/Laboratory of Environmental Chemistry and Biochemistry, Vodnany

Keywords: LDTD, psychoactive compounds, microsamples, brain, juvenile fish

Introduction

Psychoactive pharmaceuticals exhibit continuously increased consumption trends[1]. Although a concentration reached only residual levels a negative effect on aquatic organism have been observed[2]. Laser diode thermal desorption (LDTD) is direct sample introduction technique without previous analyte separation. The new analytical method with this instrumentation have been developed for psychoactive compound analysis in tissue microsamples.

Methods

The method has been evaluated for low (LDTD-QqQ) and high-resolution (LDTD-HRPS) mass spectrometers then compared with previously developed method employed liquid chromatography (LC-HRPS). The exclusion of separation technique brings some benefits such as dramatic reduction of analysis time requirement and exclusion of sample filtration as pre-treatment step which can be beneficial specially when total sample amount is limited.

Results

The LDTD parameters were carefully evaluated to achieve sufficient number of data points, the highest and the most stable response.

Linearity of calibration curve (0.1-50 ng·g⁻¹) were confirmed for all mass spectrometers and compounds by R² coefficient (no lower than 0.9950).

Developed method was validated to brain tissue samples fortified at four concentration levels. For example, average recovery obtained for citalopram at concentration of 5 ng·g⁻¹ was 102±4 % (LDTD-HRPS) and 136±9 % (LDTD-QqQ), respectively. Similar results were obtained for tramadol in a case of QqQ.

Conclusions

Validated method was used for analysis of brain tissue samples of juvenile fish exposed to psychoactive compounds at environmentally relevant concentrations. Preliminary obtained correspond with results from LC-HRPS and show potential of LDTD-QqQ and LDTD-HRPS method for analysis of complex biological samples in very short run time.

Novel Aspect

The newly developed analytical method is suitable for biological microsamples analysis for environmental research purpose.

References

1. Gualano, Maria Rosaria et al. *Psychiatric Services*, 65(10): 1226–1231 (2014).
2. Brodin, T., Fick, J., Jonsson, M., Klaminder, J., *Science*, 339:339:814-815 (2013).

Acknowledgments

CENAKVA (No. CZ.1.05/2.1.00/01.0024) and CENAKVA II (No. LO1205 under the NPU I program), Czech Science Foundation (No. 16-06498S).

MP-212 / BREATH ANALYSIS WITH SECONDARY IONIZATION HIGH-RESOLUTION MASS SPECTROMETRY (SESI-HRMS) IN THE CLINICAL ENVIRONMENT WITH A FOCUS ON DISEASE SPECIFIC BIOMARKERS IN CHILDREN WITH CYSTIC FIBROSIS

Tobias Bruderer⁽¹⁾ - ***Naemi Haas-Baumann***⁽¹⁾ - ***Ronja Weber***⁽¹⁾ - ***Astghik Baghdasaryan***⁽²⁾ - ***Srdjan Johnny Micic***⁽¹⁾ - ***Simona Mueller***⁽³⁾ - ***Florian Singer***⁽¹⁾ - ***Renato Zenobi***⁽³⁾ - ***Alexander Moeller***⁽¹⁾

University Children's Hospital Zurich, Division of Respiratory Medicine, Zurich (1) - Joint Medical Center Arabkir, Division of Allergology, Yerevan (2) - ETH Zurich, Department of Chemistry and Applied Bioscience, Zurich (3)

Keywords: high-resolution mass spectrometry, ambient mass spectrometry, clinics, children, cystic fibrosis

Introduction

This is the first study to assess real-time "breath-prints" of children with SESI-HRMS (1, 2, 3) in the clinical environment. In the context of cystic fibrosis (CF), major contributors to morbidity are early and often subclinical pulmonary infection and pronounced neutrophilic inflammation. Early detection of disease and disease associated complications is crucial for implementing timely therapeutic measures to reduce disease burden and improve prognosis.

Methods

A novel SESI ion source was used with highly inert surfaces. We present data of 41 children (23 CF and 17 healthy controls, HC), 4 - 16 years old, from an ongoing cross-sectional study. Carefully pre-processed MS data was used for biomarker detection and classification. Mann-Whitney U tests together with q-values were applied to isolate m/z values. Subsequent variable reduction through principal component analysis was used to perform binary logistic regression.

Results

Breath analysis with SESI-MS in the clinical environment requires careful control of temperatures, gas flows, contaminations (*e.g.* fatty acids from skin), cleaning procedures (*e.g.* cluster formation) and polarity switching. CF and HC had comparable BMI (17.7 ± 3.0 vs. 16.9 ± 2.3 , $p = 0.17$) and FEV1 (z-scores -0.45 ± 1.18 vs. -0.27 ± 1.13 , $p = 0.37$). Four children were *Pseudomonas aeruginosa* positive. We were able to isolate 43 m/z peaks with elevated intensity signals in CF vs. HC below an FDR control level of $q = 0.15$. Classification was applied to the first three principal components and resulted in significant outcomes (CF vs. HC) with an average accuracy of 81.3 % and a Cohen's Kappa value of 0.61. We could annotate the detected m/z peaks with molecular formulae and identify several homologues series of compounds. Structure elucidation is on-going based on comparison of fragment spectra with reference compounds for unequivocal confirmation of the involved molecular markers.

Conclusions

Non-invasive breath analysis with SESI-HRMS identifies CF specific compounds in children. This will allow further phenotyping of CF lung disease, with the potential for early detection of airway infection.

Novel Aspect

This is the first study to assess real-time "breath-prints" of children with SESI-MS in the clinical environment. We could identify CF specific volatile compounds in children.

References

1. Sinues P, Tarokh L, Li X, Kohler M, Brown SA, Zenobi R, Dallmann R. PLoS ONE, 9, 12 (2014).
2. Schwarz EI, Sinues P, Bregy L, Gaisl T, Diego-Garcia G, Gaugg MT, Suter Y, Stebler N, Nussbaumer-Ochsner Y, Bloch KE, Stradling JR, Zenobi R, Kohler M, Thorax, 71, 110-117 (2016).
3. Gaugg MT, Bruderer T, Nowak N, Eiffert L, Sinues P, Kohler M, Zenobi R, Analytical Chemistry, 89, 10329-10334 (2017).

MP-213 / ONE SIZE DOES NOT FIT ALL: DART AND COMPLEMENTARY AMBIENT IONIZATION METHODS

Robert B. Cody

JEOL USA, Inc., Mass Spectrometry, Peabody, MA

Keywords: ambient, ionization, DART, paper spray, flame

Introduction

Each of the ambient ion sources has specific strengths and weaknesses. No one method is universal. The application of several ambient ionization sources (DART, paper spray, inlet ionization and thermal ionization) to widely varying analytical problems will be presented from the perspective of the analytical chemist where no two sample types are the same.

Methods

Samples were analyzed using a high-resolution time-of-flight mass spectrometer with an exposed atmospheric pressure interface. A Direct Analysis in Real Time (DART) ion source was coupled with a thermal ionization and pyrolysis device. Ambient ionization included a prototype paper spray source, a butane torch, and inlet ionization with 3-nitrobenzotrile matrix.

Results

A surprising application is the use of DART analysis for trace quantitation of D₂O in water. DART ionization with the thermal desorption source was applied to material analysis such as forensic identification of automotive paint chips and failure analysis, yet direct analysis of natural rubber gave better results than thermal desorption. The latter observation supports the hypothesis that the DART mechanism involves more than just thermal desorption. Other applications could not be analyzed by DART. Thermal ionization was used to detect trace iodide in milk. Inlet ionization was found to be the best method for rapid quantitation of polyphenols in tea infusions, while paper spray was the best method for copolymer analysis and detection of anthocyanins in berry extracts.

Conclusions

Although the DART mechanism involves ionization of atmospheric water, DART can still be used for trace D₂O quantitation. Temperature-programmed thermal desorption/DART is useful for material analysis, but it is not always the best approach. The DART mechanism may involve more than simple thermal desorption. While paper spray and inlet ionization produce similar ion types and similar mass spectra, their applications can be complementary.

Novel Aspect

Contrasting different ambient ionization methods for unusual applications and observing some mechanistic clues

MP-214 / ENHANCED COUPLING OF ACOUSTIC LEVITATION SAIL MASS SPECTROMETRY FOR DIRECT IONIZATION AND REACTION MONITORING

Elizabeth A. Crawford⁽¹⁾ - *Demian Dietrich*⁽¹⁾ - *Cemal Esen*⁽²⁾ - *Dietrich A. Volmer*⁽³⁾

Saarland University, Institute of Bioanalytical Chemistry, Saarbruecken (1) - Ruhr-University Bochum, Department of Mechanical Engineering, Bochum (2) - Humboldt University Berlin, Department of Chemistry, Berlin (3)

Keywords: Containerless reactions, ambient mass spectrometry, acoustic levitation, reaction miniaturization and acceleration

Introduction

Acoustic (ultrasonic) levitation has emerged as the most practical method of object levitation. In our previous work using DART-MS [1], the heated gas beam of the DART ionization source could not penetrate the acoustic field of the levitator. In this work, solvent-assisted inlet ionization (SAIL), namely ambient ionization occurring directly in the heated inlet of the mass spectrometer, permitted collection and ionization of samples in a single step.

Methods

Solvent-assisted inlet ionization (SAIL) was coupled to an ion trap mass spectrometer. Online ionization of microdroplet reactions suspended under ambient conditions in an acoustic levitator was achieved directly at the inlet of the mass spectrometer. Fused silica capillary (ID 100 μm x OD 360 μm) dimensions employed in SAIL were critical for optimal ionization and external resistive heating of the capillary inlet enhanced analyte signal intensities.

Results

Solvent-assisted inlet ionization (SAIL) greatly improved the direct sampling capabilities for the coupling of liquid samples suspended in an acoustic levitator with online mass spectrometric detection. The fused silica tubing employed in SAIL proved optimal for directly extracting droplets from the levitator, simultaneously acting as the sample transfer tube and ionization source. Positioning of the fused silica capillary within the inlet also played a key role and the optimized ionization conditions were achieved by placing the ionizing end of the SAIL capillary less than a millimeter into the inlet of the mass spectrometer. A derivatization reaction of cholest-4-en-3-one with Gerard's Reagent T consisting of a 6 μL total reaction volume was compared with bulk (600 μL) reactions. Reactions carried out in the reduced reaction volume (in stillo) were approx. 10-fold faster. It was shown that external heating of the inlet of the mass spectrometer provided a 10-fold increase in signal intensity.

Conclusions

Miniaturizing the reaction scheme to the microliter scale reduced the materials needed to carry out the reaction. Performing reactions in an acoustic field allowed all of the materials to react without interferences from container walls, therefore also enhancing reaction miniaturization. Coupling the micro-reactions with SAIL-MS enabled online real time mass spectrometric monitoring.

Novel Aspect

Coupling of acoustically levitated droplets (in stillo) with SAIL mass spectrometry enabled enhanced online container less micro-reaction monitoring.

References

1. Crawford E. A., Esen C., Volmer D. A., Anal. Chem., 88, 8396-8403 (2016).

MP-215 / DIRECT INFUSION RESONANCE-ENHANCED MULTIPHOTON IONIZATION MASS SPECTROMETRY OF LIQUID SAMPLES UNDER VACUUM CONDITIONS

Hendryk Czech ⁽¹⁾ - **Claudia Kruth** ⁽¹⁾ - **Martin Sklorz** ⁽¹⁾ - **Johannes Passig** ⁽¹⁾ - **Sven Ehlert** ⁽²⁾ - **Achille Cappiello** ⁽³⁾ - **Ralf Zimmermann** ⁽¹⁾

University of Rostock, Chair of Analytical Chemistry, Rostock (1) - Photonion GmbH, Research & Development, Schwerin (2) - University of Urbino, LC-MS Laboratory, Urbino (3)

Keywords: mass spectrometry, photoionization, atmospheric pressure, petroleum, PAH

Introduction

The generation of ions with ambient ionization techniques are affected by the sample matrix leading to concurrent ion yields of molecular ions, quasi-molecular ions formed by adducts and fragments. In this study, we present direct-infusion resonance-enhanced multi-photon ionization (DI-REMPI) at 266 nm combined with a direct liquid interface [1] for soft and matrix-independent analysis of aromatics in complex oil samples.

Methods

The instrumental setup basically consists of a nano-flow UPLC for sample injection, a quadrupole mass spectrometer [2] and an Nd:YAG laser. A polycyclic aromatic hydrocarbon (PAH) solution, five diluted crude oils and four diesel samples were injected by a manual injection valve of 20 nl in a continuous solvent flow (400 nl/min) of acetonitrile. Only the tip of the capillary is heated to 280°C and ends 1 mm in front of the repeller inside the ion source.

Results

The total eluent enters the ion source where it becomes promptly nebulized followed by desolvation, which is supported by the heater enclosing the repeller. REMPI under vacuum conditions and moderate laser intensities refers to a soft, selective and sensitive ionization technique for PAHs [3], so molecular ions of six PAHs were detected from the standard solution with good linear relation between intensity and concentration. Moreover, the addition of 1, 3, 10 and 30% toluene to one crude oil did not affect the mass spectrum except its molecular ion.

In the next step, crude oils were analyzed and compared to REMPI mass spectra of evolved gases released from a thermal balance [4]. Due to the total transport of sample into the ion source by DI-REMPI, the m/z range is extended that compounds which decompose in the thermal balance can still be detected. However, homologue series of odd m/z between 200 and 300 dominated the mass spectra of the crude oils which are assigned to ionization fragments of species with high molecular weight and low ionization energies.

Conclusions

DI-REMPI with direct liquid interface enables to enhance the m/z range of analytes by introducing the entire sample into the ion source while maintaining all advantages of REMPI under vacuum conditions, including no influence of matrix, selectivity for aromatics and the generation of molecular ions. Therefore, DI-REMPI is an alternative ionization technique for liquid chromatography or direct mass spectrometry.

Novel Aspect

DI-REMPI is a novel soft ionization technique for liquid chromatography or direct mass spectrometry to investigate aromatics in liquid samples without matrix effects.

References

1. Cappiello A., Famigliani G., Pierini E., Palma P., Trufelli H., Anal. Chem. 79, 5364 (2007).
2. Schepler C., Sklorz, M., Passig, J., Famigliani, G., Cappiello A., Zimmermann R., Anal. Bioanal. Chem., 405, 6953 (2013)
3. Boesl U., J. Mass Spectrom., 35, 289 (2000)

4. Streibel T., Mitschke S., Adam T., Weh J., Zimmermann R., J. Anal. Appl. Pyrol., 79, 24 (2007)

MP-216 / ION CHEMISTRY FUNDAMENTALS OF SECONDARY ELECTROSPRAY IONIZATION MASS SPECTROMETRY (SESI-MS) FOR VOC ANALYSES

Kseniya Dryahina

J. Heyrovský Institute of Physical Chemistry, Department of Chemistry of Ions in Gaseous Phase, Prague

Keywords: SESI-MS, VOC, chemical ionization, ion-molecular reactions, SIFT-MS

Introduction

SESI-MS is a highly sensitive technique for the analyses of trace gaseous analytes in air and breath. It has been previously applied to several metabolomic studies, including lung infection in a murine model and identification of volatile metabolites with molecular weights up to 1000 Da in human breath. However, a major weakness of this method is a lack of accurate quantification of gaseous analytes, which hinders its application in the biological and medical areas.

Methods

The physical and chemical processes, which lead to gas-phase ionization of VOC in electrospray and influence of transfer line on the ion composition, are studied. Sequences of ion-molecular reactions involved in the formation of product ions are described and the chemical kinetics is assessed for equilibrium and non-equilibrium situations. Chemical kinetics modeling is exploited combining analytical and numerical methods based on experimental data.

Results

Chemical kinetics of hydrate ion formation is described as a function of temperature and partial pressure of water vapor in the gas present in the ionization chamber. The influence of transfer from the atmospheric pressure ionization region to the vacuum of the mass spectrometer via heated capillary on the ion composition is characterized. Representative mass spectral data are shown combining the SESI-MS and the selected ion flow tube mass spectrometry, SIFT-MS, experiments. Factors influencing the quantitative sensitivity of SESI-MS are summarized.

Conclusions

The step-by-step understanding of chemical and physical processes affecting the composition of ions on their journey from the ESI droplets, via the atmospheric pressure gas containing traces of analyte, through the ion transfer line into the low pressure regions of the mass spectrometer will lead into robust methods of quantification of VOCs from the SESI-MS data.

Novel Aspect

Description of sequence of processes involved in the formation of product ions in SESI will lead into a practical method of absolute quantification of trace gaseous analytes.

References

1. Zhu J. J., Jimenez-Diaz J., Bean H. D., et al. Journal of Breath Research 7, 016003 (2013)
2. Garcia-Gomez D., Sinues P. M. L., Barrios-Collado C. et al. Analytical Chemistry 87, 3087-93 (2015)
3. Singh K. D., Vidal-De-Miguel G., Gaugg M. T. et al. Journal of Breath Research 12, 027113 (2018)

MP-217 / SPIDERMASS A NOVEL TOOL TO TACKLE IN-VIVO MASS SPECTROMETRY

Philippe Saudemont⁽¹⁾ - ***Emmanuel Bouchaert***⁽²⁾ - ***Quentin Pascal***⁽²⁾ - ***Michael Ziskind***⁽³⁾ - ***Dominique Tierny***⁽²⁾ - ***Michel Salzet***⁽¹⁾ - ***Isabelle Fournier***⁽¹⁾

PRISM Inserm U1192, University of Lille, Lille (1) - Oncovet Clinical Research, University of Lille, Lille (2) - CNRS UMR 8523 PhLAM, University of Lille, Lille (3)

Keywords: in-vivo, real-time, Laser, diagnostic, surgery

Introduction

Over the past 20 years, Mass Spectrometry has been translating progressively from analyzing extracts to raw surfaces. This went with the developments of technologies operating under ambient conditions all grouped within the field of Ambient Ionization Mass Spectrometry (AIMS) [1]. AIMS enable the analysis of raw material without sample preparation but are not suited for in-vivo analysis and at least not with minimal invasiveness. We have developed a novel system for in-vivo MS analysis and in-man use.

Methods

SpiderMass [2] is composed of an OPO laser system tuned to 2.94 μm to access the most intense absorption band of water molecules. The laser irradiation promotes the excitation of water present within the tissues leading to a phase explosion. The plume of desorbed material is transferred to the MS instrument by aspiration through a tubing of several meters length connected to the MS instrument through a dedicated interface which replaces the conventional source.

Results

We have studied and optimized the performances of our instrument for different applications. In particular, we have developed a fibered laser terminated by a handpiece allowing free scanning of all surface types. At the laser energy requested for ion production, this system is of low invasiveness only leaving a temporary white trace, which disappears within minutes after the experiment, at the surface of scan tissues. We have then assessed the system for different applications. In oncology, we demonstrate that SpiderMass molecular signatures have the requested specificity and sensitivity for classification purposes of cancer types and grades as exemplified with sarcoma. Further experiments performed in-vivo on model mice developing breast cancer and then intraoperatively at the veterinary surgery room confirm the performances of SpiderMass for cancer research. But SpiderMass reveals also to be well suited for fast diagnostic from FFPE tissues at the pathologist bench, bacteria typing or dermatology with real-time in-man DMPK.

Conclusions

The novel developed SpiderMass system with its fibered laser demonstrates to be an easy and robust tool for MS analysis in-vivo with a minimal invasiveness. This system opens the door to many applications in various fields including clinics for surgeons and pathologists. The capacities of the system to operate in both positive and negative ion modes and the possible analysis of proteins makes it as well as very interesting tool for biological research.

Novel Aspect

Novel system enabling in-vivo MS analysis

References

- 1 Clendinen C. S., Monge M. E., Fernández F. M., *Analyst* 142(17), 3101-3117 (2017)
- 2 Fatou B., Saudemont P., Leblanc E., Vinatier D., Mesdag V., Wisztorski M., Focsa C., Salzet M., Ziskind M., Fournier I. *Sci. Rep.* 6, 25919 (2016)

MP-218 / MOLECULAR SURGERY NAVIGATION BY THE COMBINATION OF MASS SPECTROMETRY WITH MOLECULAR RESONANCE PROBE DESORPTION

Vladimir Frankevich - Vitaliy Chagovets - Alexey Kononikhin - Gennady Sukhikh

Federal State Budget Institution "Research Center for Obstetrics, Gynecology and Perinatology" Ministry of Healthcare and Social Development of the R, Systems Biology, Moscow

Keywords: Surgery navigation, Cancer, Molecular resonance, EESI

Introduction

Quantum Molecular Resonance (QMR) technology was recently developed as a surgical tool for tissue dissection [1]. In QMR a well-defined spectrum of high-frequency electromagnetic field is produced to selectively break molecular bonds without inducing a significant rise in local temperature. In general purposes QMR was designed for the surgery. Here we demonstrate QMR technology in combination with MS as a novel Surgery Navigation technique. We show that this technique allows on-line ambient molecular analysis of biological tissues with high chemical sensitivity.

Methods

The quantum molecular resonance generator was used as a desorption probe for breast cancer tissues chemical scanning. The probe was coupled with air jet pump which delivered desorbed molecules into the home-made EESI source coupled with qTOF mass spectrometer (Maxis, Bruker Daltonics, Germany). Tissue samples were obtained in the Research Center for Obstetrics, Gynecology and Perinatology and stored under 193 K until measurements. All data were processed with PCA and OPLS-DA methods to determine the most significant differences between tissues.

Results

The method utilized in the study allowed obtaining highly intensive, stable and reproducible signal from tissues in both positive and negative ion modes. Molecular composition differences in real time were shown for different cancer tissues. Multivariate analysis was used to classify the mass spectrometric data. Score plots from OPLS-DA revealed good clusterization of data points with accordance to tissue source. Real-time feedback to the surgeon about margin status in patients with breast cancer was shown. High accuracy tumor boundary detection is also presented.

Conclusions

A novel on-line ambient platform for biological tissue analysis is introduced

Novel Aspect

Molecular Surgery Navigation system based on the combination of mass spectrometry with molecular resonance probe desorption was developed.

References

1. Pozzato G, Vignato G: Teoria della risonanza quantica molecolare nella realizzazione del bisturi elettronico 'Vesalius'. *Quintessence Int*;5/6, 153–155, 2003
2. Chen, H., A. Venter, and R. G. Cooks, Extractive electrospray ionization for direct analysis of undiluted urine, milk and other complex mixtures without sample preparation. *Chemical Communications*, 19, 2042-2044, 2006.

This work was supported by Ministry of Education and Science of the Russian Federation grant (agreement No. 14.613.21.0059, RFMEFI61316X0059).

MP-219 / A COMPARISON BETWEEN DIATHERMIC KNIFE AND CO₂-LASER COUPLED TO REIMS-TOF FOR INSTANTANEOUS TISSUE CLASSIFICATION

Michele Genangeli - Ron M. A. Heeren - Tiffany Porta

M4I, Maastricht University, Maastricht

Keywords: , CO₂-Laser, MSI, DESI, REIMS, TOF

Introduction

Lately, numerous studies have shown the benefit of coupling rapid evaporative ionization mass spectrometry (REIMS) to a diathermic knife for rapid tissue sampling and characterization. In addition, CO₂ laser can also be used for the same purpose but in a less invasive manner. Herein, we assess the opportunities and limitations of the laser and diathermic knife coupled with REIMS based on quality of data generated and invasiveness of both sampling techniques.

Methods

Fresh tissue (*i.e.* muscle, liver, bone, bone marrow, cartilage, skin, fat) from different animals were analysed using the two sampling tools. Generated vapors were directly analyzed by a REIMS-TOF system, and tissue/species-specific molecular profiles were recorded. These profiles were integrated in a database and multivariate statistical analysis was employed to build a classification/prediction model which would function as real-time recognition software.

Results

It was possible to generate tissue vapors rich in lipids with both sampling tools. Each tissue-specific mass spectrum (= molecular fingerprint) can be recorded within 3 seconds after ablation and cauterization. The laser system allows controlling penetration power and wave type, which results in cleaner spectra with better overall signal, higher S/N ratio and higher reproducibility when compared to the diathermic knife. A statistical model employing principal component analysis and linear discriminant analysis (PCA-LDA) was built based on the characterization of different molecular profiles such as tissue type and animal species. When using one or the other technique on the same tissues, similar molecular profiles are displayed, which demonstrates the possibility to build a common database from both sampling tools.

After analysis, samples were frozen and cut into 10 µm sections. Further investigations via DESI-TOF imaging were conducted to the different lipids composition and distribution in the tissues around the burnt area.

Conclusions

Combining REIMS with different sampling interface has shown promises for several applications such as food analysis. We demonstrate laser-REIMS to produce clean cut and more reproducible signal, and the possibility to sample a wider variety of samples including hard tissues, which is not possible with the diathermic knife. The possibility to build a common database for both tools renders REIMS more versatile and extends its field of possible applications.

Novel Aspect

Novel versatile and robust applications of REIMS technology with the possibility to build a common database between sampling techniques for instant sample recognition.

References

1. Balog J. et al. Identification of the Species of Origin for Meat Products by Rapid Evaporative Ionization Mass Spectrometry, *J. Agric. Food Chem.* 64, 4793-800 (2016). doi: 10.1021/acs.jafc.6b01041.
2. Black C, Chevallier OP, et al. A real time metabolomic profiling approach to detecting fish fraud using rapid evaporative ionisation mass spectrometry, *Metabolomics* 13, 153 (2017). doi: 10.1007/s11306-017-1291-y

MP-220 / MULTI-SCHEME CHEMICAL IONIZATION INLET FOR ATMOSPHERIC PRESSURE MASS SPECTROMETRY APPLICATIONS

Matti P. Rissanen⁽¹⁾ - **Jyri Mikkilä**⁽²⁾ - **Jani Hakala**⁽²⁾

Keywords: Atmospheric pressure chemical ionization, Multiple reagents

Introduction

Atmospheric Pressure Chemical Ionization (APCI) is a sensitive and selective method to ionize the sample for mass spectrometric analysis. Different ionization chemistries ionize different compounds, and using several ionization schemes will give a comprehensive view on the composition of the sample, without losing the selectivity of the chemical ionization method.

Methods

We introduce a multi-scheme chemical ionization inlet (MION) capable of combining multiple ionization schemes in one instrument. The MION consists one or more ion sources, from which the reagent ions are mixed in the sample without contaminating the sample with the neutral reagents. The ion sources can also be turned off to measure the naturally occurring ambient ions. The MION can also be equipped with a filter desorber/calibrator unit allowing the analysis of filter samples and liquid microinjections.

Results

We studied the performance of the MION with two reagent ions, nitrate (NO_3^-) and bromide (Br^-). The MION is capable of switching between the ionization schemes within seconds, and turning the ionization sources off and switching into the ambient ion mode is as fast. The neutral reagents are prevented from mixing with the sample, leaving the sample as unperturbed as possible.

We determined the sensitivity of the MION with gaseous sulfuric acid in NO_3^- mode. Sulfuric acid represents the compounds that get ionized at (or close to) the collisional limit. The sensitivity of the MION was found to be on par with the current state of the art NO_3^- -chemical ionization inlet at the low ppq/s level [1].

We also explored the alpha pinene ozonolysis. The reaction results in, among others, a multitude of compounds having a formula $\text{C}_{10}\text{H}_x\text{O}_y$, where the amount of oxygen (y) varies over a wide range [2]. The Br^- chemistry is selective towards the less oxidized ($y = 2 - 10$) species, including the hydroperoxyl radical, whereas the NO_3^- chemistry is towards the more oxidized species ($y = 4 - 11$).

Conclusions

Our new MION inlet greatly expands the possibilities what can be done in the field of APCI-mass spectrometry with one mass spectrometer. The fast switching and superb sensitivity allows for as many ionization schemes as the user finds practical, without losing much of the time resolution.

Novel Aspect

We have produced a new atmospheric pressure chemical ionization inlet, MION, capable of switching from one ionization scheme to another within seconds.

References

1. Jokinen T. et al, Atmos. Chem. Phys., 12, 4117-4125 (2012).
2. Rissanen M. et al, J. Phys. Chem. A, 119, 4633-4650 (2015).

MP-221 / OBSERVABLE CHARGE STATES IN NATIVE MS: DIFFERENCES BETWEEN ESI AND LILBID

Jan Hoffmann - Kudratullah Karimi - Nina Morgner

University, Goethe University Frankfurt / Institute of Physical and Theoretical Chemistry, Frankfurt

Keywords: native MS, observable charge states, nESI, LILBID

Introduction

Many applications used in mass spectrometers are highly charge state dependent processes. This includes CID (collision induced dissociation), CIU (collision induced unfolding) as well as IM-MS (ion mobility MS). Different ionization methods result in different detectable charge state distributions. Depending on the application highly and lowly charged species both have their advantages and disadvantages respectively.

Methods

In LILBID [1] droplets of the analyte solution are transferred into vacuum and irradiated by IR laser pulses. This leads to an explosive expansion of the droplet and thus release of analyte ions into gas phase. The ions are collected using an ion funnel and transferred into a modified qToF [2] (QTOF I; MS Vision) mass spectrometer. Different mechanisms of dissociation (laser power and collision cell) can be combined. For comparison a commercially available nanoESI-MS instrument (Synapt G2s; Waters) was used.

Results

To investigate the observable charge state distributions samples of both native and chemically modified[3,4] lysozyme were analyzed. Modifying lysosome was achieved by acetylation (alternatively: by propylation) of amino groups by varying excess of acetic (propionic) anhydride at basic pH value. According to the protein's isoelectric point a change in charge distribution could be observed either with varying pH or with different degree of acetylation in either ion mode. Typically the actively ionizing nanoESI results in highly charged species which is in contrast to the non-actively ionizing / desorbing LILBID ion source where mostly lowly charged species are detected. In addition, with increasing hydrophobicity due to increasing modifications signal intensity issues also occurred.

Conclusions

When analyzing the average per-species charge states a decrease depending on the degree of modification could be observed as expected. Both ion sources behave very similar in this respect. Interestingly the actual decrease in charge state is significantly lower than expected indicating more complex underlying mechanisms.

Novel Aspect

LILBID coupled to a qToF enabled the investigation of small modifications to biomolecules used here to study the observable charge states of different ion sources.

References

1. Morgner N., Barth H. D., Brutschy B., Australian Journal of Chemistry, 59 (2), 109-114 (2006).
2. Sobott F., Hernandez H., McCammon M. G., Tito M. A., Robinson C. V., Analytical Chemistry, 74, 1402-1407 (2002).
3. Carbeck J. D., Severs J. C., Gao J. M., Journal of Physical Chemistry B, 102, 10596-10601 (1998).
4. Krusemark C. J., Frey B. L., Belshaw P. J., Journal of the American Society for Mass Spectrometry, 20, 1617-1625 (2009).

MP-222 / PERFORMANCE ENHANCEMENT USING HE BUFFER GAS IN A SMALL-SCALE LIT FOR IN-THE-FILED MASS SPECTROMETRY

Aurika Janulyte - Yves Zerega - Jacques Andre

Aix-Marseille University, CNRS, Laboratoires PIIM (UMR 7345), Madirel (UMR 7246), Marseille, France

Keywords: Instrumentation; in-the-field and miniature MS; LIT; Ion/buffer gas collision

Introduction

A LIT, composed of hyperbolic electrodes and disk endcaps, operates a Single Ion Monitoring mode

consisting in a series of injection, cooling, isolation and ejection steps, realized by means of only DC switched voltages, to simplify conception and guarantee stability [1,2].

The collisional process of ions with buffer gas is simulated with the aim of improving mass analysis performances depending on: peak width, analysis duration and relative number of ions involved.

Methods

The design of the electrodes, then ion trajectories have been computed using CPO software[3]. An additional subroutine has been implemented for modeling ion/buffer gas colliding as follows: random drawing according to Poisson law to define collision times (event-driven simulation); modifying ion velocity according to elastic collision law with random uniform scattering angle after collision.

Results

Collisional process does not induce a significant modification of ion state(distributions in positions and velocities) due to short duration-time of trap filling with ions at injection, with He pressures less than 10^{-3} torr.

For cooling stage, the trap is operated in RF only mode. The collisional process improves the relative number of confined ions. The cooling duration-time, to reach equilibrium ion state, has been estimated at 4.5 ms, involving almost 70 collisions, with He pressure equal to 10^{-3} torr. At equilibrium ion state, the distributions tend to Gaussian distributions, except for the distribution of axial ion positions.

For ion isolation stage, the trap is operated close to the apex. In these conditions, the cooling process induces fast-ion loss. A set of values has been found: 150 μ s for isolation step time-duration and 50% of injected ions involved in the mass spectra, with mass separation of 1u and He pressure equal to 10^{-3} torr.

Conclusions

To reduce ion loss observed during isolation step, the He buffer gas pressure could be reduced, however the duration of cooling time will be increase by the same factor. Dividing by two the pressure, the number of confined ions during isolation step has been found up to 70 % and the duration of the sequence is about 10 ms.

Mass separation could be improved by a specific ramp to move the trap operating point from the center toward the apex of the stability diagram.

Novel Aspect

Same performances have been obtained with a small-scale LIT than with a mass filter with larger length for in-the-field application.

References

1. Janulyte A.,Zerega Y.,Andre J.,Brkic B.,Taylor S., Rapid Communications in Mass Spectrometry, 30, 2407-2415 (2016).
2. Brkic B., Giannoukos S., France N., Murcott R., Siviero F., Taylor S., International Journal of Mass Spectrometry, 369,30-35 (2014).
3. Harting E., Read F. H.,Electrostatic lenses, by Elsevier Publishing Company, Amsterdam (1976).

MP-223 / CHARACTERIZATION OF COLD ELECTRON IONIZATION FOR PORTABLE MASS SPECTROMETERS

Seung Yong Kim⁽¹⁾ - *Wanseop Jeong*⁽²⁾ - *Bongyoon Yi*⁽²⁾ - *Jae Yeong Eo*⁽¹⁾ - *Byeongwon Kang*⁽³⁾ - *Hyun Sik Kim*⁽¹⁾ - *Mo Yang*⁽¹⁾

Korea Basic Science Institute, Mass Spectrometry & Advanced Instrumentation, Cheongju (1) - Korea Basic Science Institute, Mass Spectrometry & Advanced Instrumentation, Chungbuk National University, Department of Physics, Cheongju (2) - Chungbuk National University, Department of Physics, Cheongju (3)

Keywords: Cold Electron Ionization Source, Microchannel Plate, UV LED, Portable Mass Spectrometer

Introduction

Cold electron ionization (CEI) has been developed using a microchannel plate (MCP) electron multiplier initiated by UV photons from light-emitting diode (LED) [1]. Pulsed mode CEI consumed only 5 μW in average and can be controlled in nanosecond scales and is expected to replace the hot cathode electron ionization in portable mass spectrometers with the advantages such as low power consumption, low temperature, and precise control of pulse width and electron density.

Methods

A microchannel plate is used to generate the secondary electrons by illumination of an UV photons generated by various LEDs (240, 260, 280, and 300 nm wavelength) and then sequential current amplification of the electrons was performed with various combinations of electron multipliers such as microchannel plates and channeltrons. A high voltage power supply has delivered up to 2400 V. The emitted electron beam current signal is picked up by an impedance-matched anode.

Results

The property of the electron beam pulse emitted from the MCP illuminated by photon pulse from an UV LED was studied. With UV LED operated with 10 V input in 3 ns, extremely short electron beam pulse was generated which can be used as an ionization pulse in TOF MS. The measured electron beam current was more than ~ 107 electrons within 2 ns. Zubarev et al. has reported 107-108 electrons were generated by an incident nuclear fission fragments from ^{252}Cf and recorded mass spectra of various samples with a time-of-flight mass spectrometer (TOF MS) [2]. Similar ions were observed with CEI as previous reported in electron-avalanche-desorption (EAD) ionization MS. EAD ionization required sufficient electrons to desorb and ionize the molecules and record the mass spectra. For the portable TOF MS, the duration of the electron beam pulse from CEI can be precisely controlled in nanosecond scales and also provide the sufficient current, so that CEI generated enough electron to be applicable to the portable TOF MS as well as ion trap MS.

Conclusions

Cold electron ionization provided more than ~ 107 electrons within 2 ns. CEI provided high enough current to ionize the gas analyte in ion trap and TOF MS. In addition, 2 nanosecond cold electron beam pulse had been demonstrated to be an extremely short pulse for portable TOF MS. This CEIS is expected to be applicable for the portable MS with low power consumption and precise control of pulse width and electron density.

Novel Aspect

Extremely short and high current cold electron beam pulse was generated by operating a UV-LED and electron multipliers and applied for the portable ion trap and TOF MS.

References

1. H. S. Kim, S. Y. Kim, and M. Yang, *Anal Chem.*, 84, 3635-3639 (2012).
2. R. A. Zubarev, P. V. Bondarenko, A. N. Knysh, and B. V. Rozynov, *Rapid Commun. Mass Spectrom.*, 5, 278-282 (1991).

MP-224 / GRAPHENE LAYER SUBSTRATE ASSISTED LASER DESORPTION AT ATMOSPHERIC PRESSURE FOR HIGH SPATIAL RESOLUTION MASS SPECTROMETRY IMAGING

Jae Young Kim⁽¹⁾ - Sun Young Lee⁽¹⁾ - Hee Jin Lim⁽¹⁾ - Ji-Won Park⁽²⁾ - Dong-Kwon Lim⁽³⁾ - Dae Won Moon⁽¹⁾

Daegu Gyeongbuk Institute of Science and Technology (DGIST), Department of New Biology, Daegu (1) - Chungnam National University, Graduate School of Analytical Science and Technology (GRAST), Daejeon (2) - Korea University, KU-KIST Graduate School of Converging Science and Technology, Seoul (3)

Keywords: Ambient MS imaging, graphene layer, laser desorption, small molecule analysis, micro-meter spatial resolution

Introduction

We are demonstrating a new method that can facilitate the laser desorption process without gold nanoparticle treatment step. The use of graphene layer substrate is excellent platform to induce the desorption process using 532 nm-continuous wave (CW) laser. By using this method, high spatial resolution atmospheric pressure mass spectrometry (AP-MS) imaging of humid tissues can be obtained.

Methods

We used the previously developed AP-nanoPALDI MS system with a change in laser source [1]. One of the most widely used lasers, a 532 nm CW laser, was used as the desorption source. A graphene layer substrate was prepared by transferring commercial CVD graphene on Cu foil to slide glass. The specimen was placed on graphene layer substrate and positioned on the scanning stage for MS imaging.

Results

Many previous studies reported using graphene layer as a substrate for absorbing laser energy have shown that the laser source was located above the specimen and the laser light first met the specimen rather than the graphene layer[2-4]. Thus, in order for the graphene layer to role as a light absorbing layer, the thickness of the specimen had to be thin. In our proposed AP-MS system, since the laser light was focused on the specimen by the objective lens of the inverted type optical microscope, the laser source was placed under the specimen, so that the laser light was irradiated to the graphene layer first and the light energy was converted to heat from graphene layer and transferred to the tissue specimen. Therefore, micrometer-resolution ion images for mouse hippocampal tissue slices were obtained under ambient conditions.

Conclusions

Compared to the preparation of the nanoparticles treated biological specimen, while the nanoparticle treatment of the specimen was omitted, the process of transferring the graphene layer was added to the slide glass substrate. The graphene layer substrates can be prepared and stored in advance, resulting in a simplified specimen preparation and a great advantage in preparing living tissue faster.

Novel Aspect

The efficient laser desorption of bio-molecules on the wet-state hippocampal tissue is possible by use of graphene layer substrate.

References

1. Kim J. Y., Seo E. S., Kim H., Park J. -W., Lim D. -K., Moon D. W., Nat. Commun., 8, 2113 (2017).
2. Lu M., Lai Y., Chen G., Cai Z., Chem. Commun., 47, 12807-12809 (2011).
3. Lee G., Bae S. -E., Huh S., Cha S., RSC Adv. 5, 56455-56459 (2015).
4. Kim Y. -K., Na H. -K., Kwack S. -J., Ryoo S. -R., Lee Y., Hong S., Hong S., Jeong Y., Min D. -H., ACS Nano, 5, 4550-4561 (2011).

MP-225 / PORTABLE TIME-OF-FLIGHT MASS SPECTROMETER USING A COLD ELECTRON BEAM PULSE

Wanseop Jeong ^(1,2) - **Seung Yong Kim** ⁽¹⁾ - **Bongyoon Yi** ^(1,2) - **Jae Yeong Eo** ⁽¹⁾ - **Byeongwon Kang** ⁽²⁾ - **Hyun Sik Kim** ⁽¹⁾ - **Mo Yang** ⁽¹⁾

Mass Spectrometry & Advanced Instrumentation Group, Korea Basic Science Institute, Cheongju 28119,

Republic of Korea (1) - Department of Physics, Chungbuk National University, Cheongju 28644, Republic of Korea (2)

Keywords: Time-of-Flight, cold electron ionization, portable mass spectrometer, nanosecond electron pulse

Introduction

A time-of-flight mass spectrometer (TOF MS) is widely used due to fast data acquisition, wide mass range and high sensitivity. Since a very short ionization pulse is essential for high mass resolving power, a fast laser pulse is commonly used in a TOF MS. However, the size and power requirements of the laser device limit the miniaturization of a TOF MS. To solve this problem, a tiny cold electron ionization (CEI) source was suggested to generate the short ionization pulse.

Methods

The cold electrons were generated with 3 microchannel plates (MCPs) and one of MCPs was illuminated by UV photons to generate the secondary electrons [1]. They are accelerated to a metal surface where sample molecules are adsorbed. The desorbed ions are analyzed through a 18 cm ion flight tube and recorded with an oscilloscope with 2.5 GHz bandwidth. The electronic circuit for nano second electron pulse was operated with 10 V, 200 mA and 3 nanosecond.

Results

The cold electron trajectory was calculated with SIMION® Version 8.1 and analytical methods. And the calculated trajectory was confirmed with the luminescence of phosphors materials pasted on the sample plate. The cold electron ionization source recently developed in this laboratory has shown a current of up to 9 mA in 2 nanoseconds duration. This indicates the number of electrons up to 3×10^7 that would be enough for recording an usual mass spectrum. The size of cold electron ionization source is less than 2 cm^3 and the power consumption is 10 mA with a voltage of 10 V, which means 0.5 nW in average power with a 5 nanosecond per second. With the integrated CEI TOF MS, volatile gas molecules, benzene, toluene and xylene, were injected separately and detected molecular ions with H^+ and $\text{C}_n\text{H}_{2n-1}^+$ ions as previous reported in EAD (electron-avalanche-desorption) MS [2].

Conclusions

The experimental results of CEI TOF MS using a small cold electron source that generates more than 10^7 electrons in a few nanoseconds indicate that a portable CEI TOF MS with a suitable high mass resolution can detect the gas molecules. Further instrumental optimization is under development and higher electron currents and improved operating parameters will provide a useful TOF MS applicable in the field.

Novel Aspect

A portable TOF-MS was developed using an extremely short ionization electron beam pulse generated with cold electron ionization method.

References

1. H. S. Kim, S. Y. Kim, and M. Yang, *Anal Chem.* 2012, 84, 3635-3639
2. R. A. Zubarev, P. V. Bondarenko, A. N. Knysh and B. V. Rozynov, *Rapid Commun. Mass Spectrom.*, 1991, 5, 278-282

MP-226 / SOFT AND HARD IONIZATION OF VOLATILE ORGANIC COMPOUNDS BY EI FOLLOWED BY MS/MS ANALYSIS

Anna Kornilova - Lisa Cousins - Reza Javahery

Perkinelmer Inc., RND, Woodbridge

Keywords: GC-MS/MS, EI, soft and hard ionization, GC, MS, MS/MS

Introduction: Ionization Sources in Mass Spectrometry

Electron Ionization is one of the most frequently used methods in MS. While being efficient, reproducible, and non-selective, it is still a “hard ionization” method that results in high fragmentation and loss of molecular ions. A dual mode EI (DMEI) source that allows both “soft” and “hard” ionization has been developed and coupled to a GC. Such combination allows both qualitative and quantitative analyses for various volatile and semi-volatile mixtures.

Methods

Dual ion source was coupled to a modified PerkinElmer Qsight 220 triple quadrupole MS and a PerkinElmer Clarus 680 GC. Gaseous pure compounds and their mixtures were injected into the capillary column of the GC. Once separated compounds were then carried by the carrier gas (He) into the DMEI source where they were ionized; resulting ions were then analyzed by the MS/MS. System performance (reproducibility, linearity and detection limits) will be presented and discussed.

Results

Pure hexane, heptane, toluene, ethylbenzene, p-xylene, o-xylene and their mixtures were analyzed using GC-DMEI-MS/MS system. Compounds were successfully separated using GC, and then individually ionized in the DMEI. Depending on the ionization source settings two fragmentation patterns were achieved. One with minimal fragmentation and substantial molecular ion peak (“soft” ionization) and another with extensive fragmentation and almost no molecular ion peak (“hard” ionization). Different electronic settings of the ionization source allowed changing fragmentation patterns and signal intensities. Various tests showed that both modes are valuable and efficient; hence one source (DMEI) can be used for both compound identification (hard ionization) and compound quantification (soft ionization) with no change in the hardware. And since switching between hard/soft modes is almost instantaneous, both modes can be used during one sample analysis.

Conclusions

It was shown that dual mode electron ionization source can be used to perform soft and hard ionization of volatile compounds with no changes in the hardware. And while these are preliminary results and further investigations are still required, current system performance is promising and its future development is justified

Novel Aspect

introduction of an EI source that can be used in “traditional” and “soft” ionization modes. Use of GC-DMEI-MS/MS system for analysis of gaseous analytes and their mixture

MP-227 / ION SOURCE MULTIPLEXING ON A SINGLE MASS SPECTROMETER

Yury Kostyukevich - Alexander Zhrebker - Evgeny Nikolaev

Skoltech, Skoltech, Moscow

Keywords: Ion source, ESI, APPI, Isotopic exchange, Proteins, Orbitrap

Introduction

The high variety of species that need to be investigated by a mass spectrometry determined the development and use of different ion sources. For example, ESI is one of the most suitable ion sources for investigation of peptides, proteins or natural organic substances, hormones can be investigated using APCI and APPI etc. With the introduction of the ion funnel, that can effectively gather ions, it became possible to

develop a multiple inlet systems in which ions from different sources are transported into the fore vacuum region through different inlet capillaries.

Here we present the simple approach for the combination of different ion sources on a single mass spectrometer and simultaneous performing different gas phase reactions without any interference between them [1-3].

Methods

All experiments were performed on an orbitrap QExactive mass spectrometer (Thermo) with installed ion funnel and fore vacuum MALDI (Spectrograph). The front panel with the MALDI translational stage was replaced with specially developed plate with several inlet capillaries. For the experiment with multiple ion sources copper tubes were connected to inlet capillaries. Valves enabled switching channels on and off. For experiments with thermal dissociation stainless steel tubes were connected to inlet capillaries. One of the tubes was wrapped in a special heat resistant cable with a nichrome conductor in the magnesium oxide and steel braid. This cable allowed heating of the extended capillary to 600°C. To perform H/D exchange a droplet of D₂O was placed below the capillary inlet.

Results

By using separate inlet capillaries and relatively long transfer tubes with valves we have succeeded in combining 4 different ion sources: APPI, radioactive ion source, native ESI and regular ESI. Switching between channels takes less than 1s and no interference between channels was observed. We were able to observe simultaneously native tetramers of streptavidin, unfolded ubiquitin, crude oil and residue organic present in air. The developed source was used for internal calibration of humic substances by infusing clusters of phosphorus acid through another channel. Using the independence of different channels it is possible to develop a system in which ions (produced by the same or different ion sources) would be subjected to the ion-molecular reactions, such as thermal fragmentation or H/D exchange in such a way that through one channel always pass a non-modified parent ions.

We have developed a novel concept of MS-SWATH in which parent ions are observed simultaneously with the fragment ions in the same mass spectrum. Peptides of Cytochrome C digest were separated by liquid chromatography, ionized by ESI and delivered to mass spectrometer through 2 different capillaries one of which was heated to induce thermal dissociation. Another channel remains cold. By regulating flow in channels it is possible to regulate intensities of parent and fragment ions which are present in spectrum simultaneously.

Instead of thermal dissociation other gas phase reactions can be performed in different channels independently. We demonstrated this for the gas phase H/D exchange reaction. Using multiplexed ion source it is possible to rapidly switch between the channel in which H/D reaction takes place and the reference channel. It was found that channels do not affect each other.

Conclusions

We have developed a novel approach to the combination of different ion sources on a single mass spectrometer. By using separate inlet capillaries and relatively long transfer tubes with valves we have succeeded in combining 4 different ion sources: APPI, radioactive ion source, native ESI and regular ESI. Switching between channels takes less than 1s and no interference between channels was observed. Such different objects as crude oil and streptavidin tetramers can be infused and recorded simultaneously. The possibility to transport ions through a 1.5 long tube makes it possible to assemble different ion sources on different working places that would not interfere with each other.

Use of our approach considerably expands the capability of gas phase chemistry to the investigation of structure of unknown sample. Indeed, ions produced from the same ion source, possibly, preselected by FAIMS, could be split into different transfer tubes in which different gas phase reactions occur. One of the channels may be remained as reference. We have demonstrated such idea for thermal dissociation and gas phase H/D exchange reactions. Using valves it is possible to rapidly switch between different channels and it was proved that channels do not affect one another. Comparison of mass spectra recorded when different channels are open, will provide structural and chemical information about unknown species.

Novel Aspect

A multiplexed ion source allowing simultaneous combination of various ion sources and performing gas phase reactions was developed.

References

1. Y Kostyukevich, E Nikolaev Ion source multiplexing on a single mass spectrometer *Analytical chemistry* 90, 3576-3583 (2018)
2. YI Kostyukevich, AS Kononikhin, IA Popov, EN Nikolaev Analytical description of the H/D exchange kinetic of macromolecule. *Analytical chemistry*, 90, 5116–5121 (2018)
3. Y Kostyukevich, T Acter, A Zhrebker, A Ahmed, S Kim, E Nikolaev Hydrogen/deuterium exchange in mass spectrometry *Mass spectrometry reviews* (2018)

MP-228 / IDENTIFICATION OF CHEMICALS IN DIFFERENT OBJECTS BY DIRECT ANALYSIS IN REAL TIME MASS SPECTROMETRY

Ákos Kuki - Borbála Antal - Lajos Nagy - Tibor Nagy - Georgina Zelei - Miklós Zsuga - Sándor Kéki

University of Debrecen, Department of Applied Chemistry, Debrecen

Keywords: direct analysis in real time mass spectrometry (DART-MS), semi-quantitative analysis, textile articles, healthcare products, polyurethane

Introduction

The manufacture of textile, healthcare and polyurethane (PU) products employs a diverse range of processes using a large number of chemicals. Residues of these substances may remain in the products. The investigation of chemicals in these clothing and hygiene and PU articles (*e.g.* socks, shirts, sanitary articles, sterile gauze bandages, sponges, foam toys) is important since they can get into direct contact with the human skin.

Methods

Measurements were performed with a MicroTOF-Q type Qq-TOF MS instrument (Bruker Daltonics, Bremen, Germany) equipped with a DART SVP ion-source (IonSense, Saugus, MA, USA). The solid samples were manually introduced into the DART gas stream without any sample pretreatment. The results were confirmed by GC-MS and HPLC-MS.

Results

The presence of more than 100 chemicals have been revealed in 50 different commercially available clothing, personal care, healthcare, and PU items [1,2]. The identification was confirmed by DART MS/MS experiments for 42 compounds by combining the MS/MS spectra with findings from the literature and databases of spectra. In the clothing and healthcare items the most notable hazardous substances were nonylphenol ethoxylates (NPEs), phthalic acid esters, amines released by azo dyes, quinoline derivatives, tributyl phosphate, tris(2,4-di-tert-butylphenyl) phosphite, erucamide. Residues of catalysts, stabilizers, antioxidants, flame retardants, plasticizers, chain extenders, chain terminators, polyols, solvents, degradation products and contaminants, a few of them presumably toxic, were detected and identified in the PU items. Automated data acquisition and processing method and an efficient semi-quantitative DART-MS analysis were developed for the recognition, and quantitation of NPE residues, respectively.

Conclusions

The large number and diversity of the detected compounds proves unambiguously the effectiveness of DART-MS and DART-MS/MS for the rapid identification of compounds in clothing, personal care, healthcare, and PU items without any sample preparation. Our experiments might have drawn attention to

the presence of many harmful substances on the surface of commonly used products that can come into direct contact with human skin, or even blood.

Novel Aspect

Our work has explored a new application field of the DART-MS technique. A semi-quantitative DART-MS analysis was developed to determine the concentration of NPE.

References

1. Antal B., Kuki Á., Nagy L., Nagy T., Zsuga M., Kéki S., *Anal. Bioanal. Chem.*, 408, 5189-5198 (2016)
2. Kuki Á., Nagy L., Nagy T., Zsuga M., Kéki S., *Anal. Bioanal. Chem.*, 409, 6149-6162 (2017)

Acknowledgement

The work was supported by the GINOP-2.3.2-15-2016-00041 project and by the grant K-116465.

MP-229 / EXPLOSIVES DETECTION WITH A MOBILE THERMAL DESORPTION DIRECT ANALYSIS IN REAL TIME MASS SPECTROMETRY SYSTEM (TD- DART-MS)

Frederick Li - Brittany Laramee - Paul Liang - Brian Musselman

IonSense, Inc., Chemistry, Saugus

Keywords: DART-MS, Thermal Desorption, Explosives Detection, Swipe Sampling, Mobile MS System

Introduction

A thermal desorber is configured to a DART-MS system on a mobile platform to permit the detection of explosives from swipe sampling. The mobility aspect of this configuration can permit mass spectrometry based detection of explosives in the field. This work presents the capability of the system and serves as a means of determining the readiness level of the system as an explosive trace detector.

Methods

A thermal desorption unit is coupled to a DART equipped single quadrupole mass spectrometer, which is configured to a mobile cart. Standards of explosives were deposited onto PTFE-coated fiber glass swipes for analysis. A design of experiment was employed to evaluate conditions that affect detection and sensitivity. Limit of detection was determined for each explosive.

Results

Explosives were deposited directly onto PTFE-coated fiberglass swipes and analyzed using the thermal desorber operated at three temperatures: 100°C, 200°C and 300°C. Preliminary results showed the response was at the highest at 200°C for more than half of the explosives. Desorption is near instantaneous for all 10 explosives, typically requiring 5 seconds after insertion of the swipe. Complete desorption of the explosive from the swipe can be achieved in 6-12 seconds, which is dependent on the vapor pressure of the explosive and quantity on the swipe. Thermal desorber temperature, the DART grid voltage and flow rate of the pump used for DART had a significant effect on sensitivity. Sensitivities between 5 ng and 50 ng for both organic explosives such as TNT, RDX, HMTD and TATP and inorganic explosives such as potassium perchlorate were consistently achieved.

Conclusions

Detection of all explosives was consistently achieved at levels between 5 ng and 50 ng. TD-DART-MS provides a means for rapid detection of explosives from swipe sampling. The mobility aspect of this system is attractive for field detection of explosives. Future work will include determining the sensitivity of additional explosives such as smokeless powder, C-4 and homemade explosives (HMEs).

Novel Aspect

TD-DART-MS on a mobile platform can enable field detection of a wide range of explosives from swipe sampling.

MP-230 / PREDICTING AGENTS THAT “SUPERCHARGE”

Rachel Ogorzalek Loo⁽¹⁾ - Reid O'Brien-Johnson⁽²⁾ - Michael Nshanian⁽²⁾ - Joseph A. Loo⁽²⁾

UCLA, DOE Laboratory for Genomics & Proteomics, Los Angeles (1) - UCLA, Chemistry & Biochemistry, Los Angeles (2)

Keywords: supercharging, electrospray ionization mechanism

Introduction

Supercharging is the ability to increase analytes' average charge state by supplementing electrospray solvents with a low volatility compound. The supercharging methodology is useful for MS and MS/MS, but the mechanism(s) responsible for increasing charge are contentious. We describe a “supercharging” mechanism that applies to native and denaturing solutions and demonstrate its utility to predict agents that sub- and supercharge.

Methods

ESI mass spectra were acquired on hybrid quadrupole time-of-flight (QTOF; Waters Synapt G2-Si), ESI-TOF (Agilent), and 15-Tesla Fourier transform-ion cyclotron resonance (FT-ICR; Bruker Solarix) spectrometers. NanoESI sources used borosilicate glass capillaries with metal coatings (Proxeon Biosystems) operated at low analyte flow (< 50 nL/min). Final protein concentrations for ESI-MS (non-denaturing) were ~ 5 μM in 20 mM ammonium acetate, pH 6.8.

Results

Multiple mechanisms have been proposed to explain how low-volatility additives increase charge in ESI, but they typically explain only subsets of experimental conditions or analytes; *e.g.*, surface tensions exceeding those of CH₃CN or CH₃OH were credited with some additives' abilities to increase protein charge from denaturing solutions, while increases delivered by the same compounds from near-native solutions were attributed to protein denaturation, driven by the additives' increasing concentration within the evaporating droplet. The former proposal was subsequently countered by evaporation rate arguments addressing surface tension changes in evolving droplets, while the latter proposal was countered by many experimental observations.

Many factors impact charge state distributions (CSDs), but we think they can be unified into a single mechanism, applicable to a broad range of molecules from both denaturing and non-denaturing solutions. A correct mechanism will predict how new reagents affect charging in ESI.

Conclusions

The hypothesis that supercharging and sub-charging molecules have specific Bronsted acidity and basicity properties (ideally neutral molecules that are both very weak bases and very weak acids compared to H₂O), was tested by examining how adding related compounds of known pKa/pK_B (pyrazole, 1-nitropyrazole, 3-nitropyrazole, pyrrole, 1-ethyl pyrrole, 2-methyl pyrrole, and others) to ESI solvents altered CSDs. Analyte charging shifted in accordance with predictions.

Novel Aspect

Supercharging and sub-charging additives for positive ion mode can be predicted based on the solution phase basicities of low-volatility compounds.

MP-231 / SUBSTRATE MODIFICATION FOR THE ANALYSIS OF HIGHLY VOLATILE CHEMICAL WARFARE AGENTS BY PAPER SPRAY MASS SPECTROMETRY

Phillip Mach⁽¹⁾ - *Elizabeth Dhummakupt*⁽²⁾ - *Paul Demond*⁽¹⁾ - *Daniel Carmany*⁽¹⁾ - *Michael Busch*⁽¹⁾ - *Trenton Tovar*⁽²⁾ - *Ann Ploskonka*⁽³⁾ - *Jared DeCoste*⁽²⁾ - *Trevor Glaros*⁽²⁾

Excet, Inc., Research & Technology Directorate, Springfield (1) - Edgewood Chemical Biological Center, Research & Technology Directorate, Aberdeen Proving Grounds (2) - Leidos, Inc., Research & Technology Directorate, Aberdeen Proving Grounds (3)

Keywords: Paper spray, chemical warfare agents, volatility, metal-organic framework, derivatization

Introduction

Paper spray mass spectrometry (PS-MS) has been shown successfully analyze chemical warfare agent (CWA) simulants [1,2]. However, transitioning to actual G-series CWAs (*i.e.* sarin, soman) for analysis via PS-MS proves challenging due to volatility differences between the simulants and actual CWAs. In order to address this challenge, methods to fix these chemistries utilizing 2-[[dimethylamino)methyl]phenol (2-DMAMP) and metal-organic frameworks (MOFs) to extend their analytical lifetime were explored.

Methods

Two separate methods of modification were utilized: 1) MOFs were weighed out in 1:1 ratio with cornstarch and dissolved in 5 mL of water. An aliquot was spotted onto cut glass fiber filters and allowed to dry. 2) Two microliters of a 2-DMAMP solution in methanol were spotted onto paper and allowed to dry. For each method, quantities of CWAs were applied in solution or vaporized to the treated substrates. Time courses and calibration curves were subsequently carried out.

Results

Attempts to analyze G-series CWAs on paper substrates were not successful, as the agents volatilize off the paper in less than five minutes. Analysis on glass substrates only extends the analytical lifetime for up to 15 minutes. The addition of the MOFs HKUST-1 and UiO-66 to the glass substrate was found to extend the analytical lifetime of GB, GD, and GF for up to at least 60 minutes.

Paper substrates treated with 2-DMAMP were spotted with either GB or GF and analyzed over time. The G-DMAMP complexes formed were persistent out to 96 hours (the last time point tested). Additionally, calibration curves of the two G-DMAMP complexes were generated, and the limit of detection (LOD) for GB-DMAMP was 0.941 ng/mL and GF-DMAMP was 5.14 ng/mL.

Conclusions

Substrate modifications have proven useful to (1) prolong the analytical lifetime of highly volatile native G agents and (2) complex with the G agents to reduce their volatility. The MOF approach allowed for a four-fold extension in time, while preserving the ability to analyze the native molecule. The 2-DMAMP allowed for long term capture through derivatization. Sub-1 ng/mL limits of detection were achieved with G-DMAMP complexes, and the complexes remained stable for up to 96 hours.

Novel Aspect

First on-substrate modification and reaction of highly volatile chemical warfare agents.

References

1. Dhummakupt, E. S., McKenna, J., Connell, T., Demond, P. S., Miller, D. B., Nilles, J. M., Manicke, N. E., Glaros, T., *Analyst*, 142, 1442-1451 (2017).
2. Dhummakupt, E. S., Mach, P. M., Carmany, D., Demond, P. S., Moran, T. S., Connell, T., Wylie, H. S., Manicke, N. E., Nilles, J. M., Glaros, T., *Analytical Chemistry*, 89, 10866-10872 (2017).

Funding for this project is provided by the Defense Threat Reduction Agency – Joint Science and Technology Office for Chemical and Biological Defense

MP-232 / LASER-ABLATION RAPID EVAPORATIVE IONIZATION MASS SPECTROMETRY (LA-REIMS) FOR LIPIDOMICS ANALYSIS OF LIVING CHO CELLS DIRECTLY FROM THE CELL CULTURE

Stefania Maneta-Stavarakaki⁽¹⁾ - *Simon J. Cameron*⁽¹⁾ - *Julia Abda*⁽¹⁾ - *Alvaro Perdones-Montero*⁽¹⁾ - *Yuen-Ting Chim*⁽²⁾ - *Paloma Diaz-Fernandez*⁽²⁾ - *Zoltán Takáts*⁽¹⁾

Department of Surgery and Cancer, Sir Alexander Fleming Building, Imperial College, London SW7 2AZ, UK (1)- Biopharm Upstream Process Research, R&D Platform Technology & Science, GSK Medicines Research Centre, Stevenage, Hertfordshire, SG1 2NY, UK (2)

Keywords: Mass Spectrometry, Lipidomics, REIMS, Cells

Introduction

Lipidomics analysis of living cells is still a difficult task, mainly because of the low abundance lipids that demand the use of time-consuming lipid extraction methods. In the present study we report the application of LA-REIMS directly in the cell culture, in ambient environment.

Methods

LA-REIMS was performed in CHO-S cell culture liquid aliquots using the OPO (optical parametric oscillator) laser, which was connected through a REIMS source interface to a Waters G2-XS Q-ToF. The advantage of the OPO laser in contrast to others is the ability to vary the wavelength in the mid- and far-infrared region, which allows selective excitation of different functional groups.

Results

LA-REIMS connected to a Q-ToF and operated in full scan mode was able to detect most of the lipids related to CHO-S cell line, as reported in the literature [1,2], with good sensitivity.

The culture medium was used as control sample in order to distinguish the peaks coming from the cells. CHO-S cell pellets were also tested for comparison purposes. Different laser settings were tested for optimization, and the best signal was obtained with a mid-range laser power and 2900 nm wavelength.

The background “noise” was consistently low, allowing detection of low abundant lipids.

The majority of the m/z values detected in the cell culture and the cell pellet but not in the control sample were matched to putative lipids of all classes using accurate mass annotation and statistical analysis showed clear differences between the control and the sample.

Conclusions

The results of this study indicate that LA-REIMS can be used for fast and high-throughput lipidomics analysis in cell cultures without requiring any sample preparation. It is a fast and user-friendly method that is applicable in both the industry and academia for cell analysis at different levels of development.

Novel Aspect

LA-REIMS can be used for fast and high-throughput analysis of living cells as a biological snapshot of cell cultures in their native state.

References

1. Zhang Y., Baycin-Hizal D., Kumar A., Priola J., Bahri M., Heffner K., Wang M., Han X., Bowen M., Betenbaugh M., *Analytical Chemistry*, 89, 1477-1485 (2017).
2. Brugger B., Erben G., Sandhoff R., Wieland F., Lehmann W., *Cell Biology*, 94, 2339-2344 (1997).

MP-233 / OPTICAL CHARACTERISTICS OF MATRIX MATERIALS FOR MATRIX-ASSISTED LASER DESORPTION/IONIZATION

Noritaka Masaki - Shigetoshi Okazaki

Department of Medical Spectroscopy, Hamamatsu University School of Medicine, Hamamatsu

Keywords: matrix-assisted laser desorption/ionization, photo-dissociation, photosensitization, excited triplet state

Introduction

Finding new matrix materials for matrix-assisted laser desorption/ionization (MALDI) is an analytical demand for biological study to analyze wider variety of molecular species and to improve ion yield. Matrix materials have been opted conventionally based on acid–base neutralization mechanism [1]. In the present study, we aimed to characterize matrix materials from optical aspects.

Methods

Absorption spectra of matrix materials dissolved in MeOH or sprayed on quartz slide were measured with U-3500 spectrometer (Hitachi). Photo-dissociation of matrix materials was evaluated with time course of absorption spectra of matrix solutions after various irradiation periods of Nd:YAG laser. Photosensitization effect of matrix materials was detected as 1270 nm emission from singlet oxygen [2].

Results

α -Cyano-4-hydroxycinnamic acid, sinapic acid, 2,5-dihydroxybenzoic acid, 2,5-dihydroxyacetophenone, norharmane, 4-nitroaniline, 2',4',6'-trihydroxyacetophenone monohydrate, 9-aminoacridine, and 1,5-diaminonaphthalene were analyzed as typical matrices. Absorption profile obtained in a film state was more suited for wavelength at 337 nm or 355 nm than in a solution state. Absorption profile after irradiation of 355 nm Nd:YAG laser showed decrease of peaks characteristic to matrix materials reflecting photo-dissociation. As the photo-dissociation was reduced by bubbling with nitrogen gas, it dues to dissolved oxygen. In addition, photosensitized singlet oxygen was detected as 1270 nm emission from most of typical matrix materials against 355 nm Nd:YAG laser irradiation. These suggest that the laser excited matrix materials to a triplet state through intersystem crossing.

Conclusions

Photo-dissociation and/or photosensitization effect were found in all the typical matrix materials used in this study. Both photo-dissociation and photosensitization are related to production of singlet oxygen by 355 nm Nd:YAG laser irradiation. Therefore, matrix material for MALDI can be characterized as excitation to triplet state by MALDI laser or photosensitization ability.

Novel Aspect

Excitation to a triplet state with 355 nm laser was revealed as a common optical characteristic to matrix materials for MALDI.

References

1. Shroff R., et al., PNAS, 106, 10092–10096 (2009).
2. De Rosa M. C., Crutchley R. J., Coord. Chem. Rev., 233–234, 351–371 (2002).

MP-234 / EXTENDING THE POLARITY RANGE FOR DIRECT SPME-MS ANALYSES BY USING NOVEL CAPILLARY DBDI AND APPI SOURCES

Mario Francesco Mirabelli - Anna K. Huba - Renato Zenobi

ETH Zurich, Department of Chemistry and Applied Biosciences, Zurich, Switzerland

Keywords: SPME-MS, DBDI, APPI, reactivity, electrodes

Introduction

Although direct interfaces of solid-phase microextraction (SPME) with mass spectrometry (MS) are rapidly spreading and proved themselves as a fast and convenient way to analyze complex sample matrices, some limitations still need to be addressed: the limited range of compound polarity achieved by most ambient ionization sources, the sometimes severe matrix effects, and the lack of proper sample introduction strategies for a simple and fast analytical workflow

Methods

The capillary atmospheric pressure photoionization (cAPPI) and dielectric barrier discharge ionization (DBDI) sources used in this study were built in-house, in a capillary format directly connected to the MS using a leak-tight connection. For both sources, dopants were used to enhance ionization efficiency. In addition, the ionization efficiency for nonpolar compounds was increased for DBDI by using customized electrode geometry (spring- and ring-shaped)

Results

Most ambient ionization sources used for direct SPME-MS analyses are mainly sensitive towards polar compounds. Capillary DBDI and APPI allowed to efficiently ionize nonpolar compounds (*e.g.*, PAHs), introduced into the sources by thermally desorbing SPME fibers. In most cases, more polar compounds could be ionized efficiently as well. The results show a sensitivity in the pg/mL concentration range, even in the presence of relatively complex sample matrices like lake water or food samples. The analytical workflow was completely automated with a commercial PAL RTC autosampler (CTC Analytics), with a throughput of 1 sample every 1.5 minutes. For cAPPI, best results were obtained with chlorobenzene, fluorobenzene and acetone as dopants. For DBDI, the use of dopants affected the ionization mechanism, resulting in the preferential formation of radical versus protonated cations. The energy involved in the ionization processes was quantified by using substituted benzylamines thermometer ions

Conclusions

The proposed approach allowed to quantify polar and nonpolar compounds from relatively complex matrices within minutes. The range of compound polarity accessible to SPME-MS direct analyses was broadened. The selectivity of the ionization process could be tuned by using different dopants. In addition, specific DBD electrode configurations allowed to extend or contract the plasma generated inside the source, promoting chemical reactions inside the plasma

Novel Aspect

No chromatographic analysis of complex samples by efficient capillary plasma- and photoionization sources. Reactivity tuned by dopants and custom electrode configuration

MP-235 / COMBINED INFRARED AND ULTRAVIOLET LASER ABLATION ELECTROSPRAY IONIZATION MASS SPECTROMETRY

Kermit K. Murray - Remi Lawal - Fabrizio Donnarumma

Louisiana State University, Department of Chemistry, Baton Rouge, USA

Keywords: Laser, ablation, electrospray, infrared, ultraviolet

Introduction

We are using combined mid-infrared and deep-ultraviolet two-laser ablation coupled with electrospray ionization for ambient mass spectrometry of biomolecules in tissue. The goal is to increase nanoparticle

production and improve sensitivity by using the UV laser to disrupt the tissue structure followed by IR ablation and ionization.

Methods

In this work, we are using a 193 nm ArF excimer laser to disrupt the tissue prior to irradiating with a 3000 nm IR optical parametric oscillator. Both lasers are focused onto the same target spot and separated in time by an adjustable delay. A nanospray needle is directed at the inlet on-axis of a modified quadrupole time-of-flight mass spectrometer.

Results

The two-laser ablation system has been constructed and initial studies carried out for optimization of the system with peptide and protein standards. The lasers are mounted on an aluminum breadboard adjacent to the ion source and are focused onto the target with a single calcium fluoride lens for each beam. The dual-laser configuration can be operated either with the UV firing first to disrupt the covalent bonding in the tissue or with the IR firing first to heat the tissue. Initial studies with 193 nm laser ablation electrospray ionization demonstrate that the deep-UV is a much softer ionization method than might be anticipated and can produce ions from peptides and proteins without a matrix and with little fragmentation. Based on this interesting new result, initial experiments are aimed at improving the efficiency of the deep-UV ablation using IR laser pre-heating.

Conclusions

We have demonstrated that deep-UV and IR ablation coupled with electrospray ionization is a promising soft ionization method for large molecules with applications to tissue imaging. Continuing experiments are aimed at optimizing the UV and IR laser pulse energies and time delay to improve the sensitivity as well as improving the UV laser focus to improve spatial resolution.

Novel Aspect

Combined mid-infrared and deep-ultraviolet laser ablation on the same spot for laser ablation electrospray ionization imaging.

MP-236 / RAPID QUALITATIVE ANALYSIS OF TWO FLAVONOIDS, RUTIN AND SILYBIN, IN MEDICAL PILLS BY DIRECT ANALYSIS IN REAL TIME MASS SPECTROMETRY (DART-MS) COMBINED WITH IN SITU DERIVATIZATION

Lajos Nagy - Tibor Nagy - Ákos Kuki - Miklós Zsuga - Sándor Kéki

University of Debrecen, Department of Applied Chemistry, Debrecen

Keywords: direct analysis in real time mass spectrometry (DART-MS), silylation, in situ derivatization, flavonoid derivatives, silybin, rutin

Introduction

Flavonoids belong to the family of benzo-pyrone having a wide range of biological and pharmacological effects, including anti-allergic, anti-inflammatory, antioxidant, anti-cancer activities. The conventionally used mass spectrometric and chromatographic methods are quite time consuming and labor intensive. In contrast, direct analysis in real time mass spectrometry (DART-MS), enable the examination of solid samples without sample preparation.

Methods

Measurements were performed with a MicroTOF-Q type Qq-TOF MS instrument (Bruker Daltonics, Bremen, Germany) equipped with a DART SVP ion-source (IonSense, Saugus, MA, USA). Derivatization of

the standard samples was carried out on DART sampling rods by adding 2 μL of the silylation reagent to the standard powder. These samples were manually inserted into the helium gas stream.

Results

Direct analysis in real time mass spectrometry (DART-MS) with in situ silylation was used for the rapid analysis of the flavonoids (silybin, rutin) [1]. Three commercially available silylation reagents were applied. The flavonoids were detected with various degrees of silylation, and the formation of dimers was also observed. HMDS/TMCS/pyridine was the best choice for the DART-MS analysis of silybin, and BSA/TMCS/TMSI was the most effective for the detection of rutin. The effects of the DART source temperature on desorption, ionization, in-source fragmentation, dimer formation, and hydrolysis of the trimethylsilyl groups were also studied. In addition, the collision-induced dissociation (CID) properties of the derivatized silybin and rutin were explored. With our in situ silylation method, the derivatized bioactive compounds in intact medical pills could also be detected by DART-MS.

Conclusions

Our experiments showed that in situ silylation enables the use of DART-MS for the analysis of compounds that have low volatilities and/or are thermally unstable, such as silybin and rutin. Three silylating reagents were tested, and significant differences were found in their efficiencies of silylation. HMDS/TMCS/pyridine was the best choice for the DART-MS analysis of silybin, whereas BSA/TMCS/TMSI was the most effective for the detection of rutin.

Novel Aspect

Our work has explored a new application field of the DART-MS technique. A method was developed for the rapid analysis of silybin and rutin by DART-MS using in situ silylation.

References

1. Nagy T., Kuki Á., Nagy L., Zsuga M., Kéki S., J. Mass Spectrom., 53, 240-246 (2018)

Acknowledgement

The work was supported by the GINOP-2.3.2-15-2016-00041 project and by the grant K-116465 and supported through the New National Excellence Program of the Ministry of Human Capacities, ÚNKP-16-3 (T. Nagy).

MP-237 / A NOVEL MODULAR AND VERSATILE TOF PLATFORM COMBINABLE WITH VARIOUS ION SOURCES

Alfons Jordan ⁽¹⁾ - **Felix Piel** ^(1,2) - **Stefan Feil** ⁽¹⁾ - **Eugen Hartungen** ⁽¹⁾ - **Gernot Hanel** ⁽¹⁾ - **Jens Herbig** ⁽¹⁾ - **Klaus Winkler** ⁽¹⁾ - **Rene Gutmann** ⁽¹⁾ - **Lukas Märk** ⁽¹⁾ - **Philipp Sulzer** ⁽¹⁾

Ionicon Analytik GmbH, Innsbruck (1) - University of Innsbruck, Institute for Ion Physics and Applied Physics, Innsbruck (2)

Keywords: TOF, TOF-MS, mass spectrometry, APi, PTR-MS

Introduction

Time-Of-Flight (TOF) mass spectrometers are state-of-the-art analyzers not only in Proton-Transfer-Reaction – Mass Spectrometry (PTR-MS) [1], but many other fields. Depending on the application the TOF instruments range from bulky high-resolution to compact low-resolution models. Here we present a novel modular TOF design, which can be adjusted in dimensions and mass resolution and used in combination with various ion sources.

Methods

The newly developed TOF recipient consists of an orthogonal injection "base unit" equipped with an ion mirror and microchannel plate detector. The flight path of ions can be extended by installing one or two recipient extension modules, respectively. An ion funnel and a hexapole ion guide can be integrated for ion focusing, either in combination or separately. We tested the setup with a PTR source and a novel Atmospheric Pressure interface (APi) front-end.

Results

Utilizing a well established PTR ion source [2] we evaluated the performance of the modular TOF platform (0, 1 and 2 flight path extension modules, respectively) in combination with an RF ion funnel and hexapole ion guide, while analyzing a certified gas standard (TO-14A). The mass resolution of the compact base unit without any RF focusing device is above 1500 $m/\Delta m$ and the sensitivity is between 40 (below m/z 100) and more than 100 cps/ppbv (above m/z 100). By the implementation of an ion funnel the sensitivity gets significantly improved, *i.e.* to more than 500 cps/ppbv. After installation of 2 extension modules, an ion funnel and a hexapole ion guide a mass resolution of $>6000 m/\Delta m$ and sensitivities of >1000 cps/ppbv can be achieved. The latter means that after an integration time of one minute the detection limit of the instrument is below 1 pptv.

Exchanging the PTR source with an APi front-end, consisting of differential pump stages and an additional hexapole ion guide, we are able to investigate ions in the atmosphere over broad m/z range.

Conclusions

The results obtained with the modular TOF platform in combination with a PTR ion source indicate an excellent performance, which surpasses known PTR-TOFMS instruments [3]. The design is highly versatile, so that *e.g.* using two TOF base units in opposite polarity with APi front-ends, creates a revolutionary APi-TOF instrument for simultaneously detecting positive and negative ions over a broad m/z range.

Novel Aspect

A modular TOF design enables maximum flexibility utilizing extension modules while performing at high mass resolution and sensitivity.

References

1. A. M. Ellis, C. A. Mayhew, John Wiley & Sons, Chichester, UK (2014).
2. P. Sulzer, A. Edtbauer, E. Hartungen, S. Jürschik, A. Jordan, G. Hanel, S. Feil, S. Jaksch, L. Märk, T. D. Märk, *Int. J. Mass Spectrom.* 321 (2012) 66-70.
3. A. Jordan, S. Haidacher, G. Hanel, E. Hartungen, L. Märk, H. Seehauser, R. Schottkowsky, P. Sulzer, T. D. Märk, *Int. J. Mass Spectrom.* 286 (2009) 122-128.

MP-238 / IONIZATION PRODUCTS OF TERPENE SPECIES IN LTP FT-ICR-MS MEASUREMENTS

Björn Raupers⁽¹⁾ - ***Jürgen Grotemeyer***⁽¹⁾

Institute of Physical Chemistry, University of Kiel, Kiel (1)

Keywords: LTP, ICR MS, Terpenes, Cluster, Atmospheric conditions

Introduction

Previous studies [1;2] have shown signs of terpenes forming large clusters in low temperature plasma (LTP) ion mobility spectrometry hyphenations. In this work we present further studies of terpene systems using a Bruker Fourier transform ion cyclotron resonance mass spectrometer (FT-ICR-MS).

Methods

Mass analysis were performed on a 7.05 T FT-ICR system, provided by Bruker Daltonik. A LTP torch built

according to the design of Nørgaard et al. [3] with slight variations was used as an ion source. Plasma gas and sample gas were Argon (99.999 %) by Air Liquid Deutschland GmbH. As a test substances a great variety of diterpenes and sesquiterpenes were used. All purchased from Sigma Aldrich (Taufkirchen).

Results

Mass spectra of LTP FT-ICR hyphenations mostly show a great variety of ions. Predominantly Dehydrated species, oxides and fragments can be detected. Furthermore clusters of all these can be found. Depending on the analyte species clusters can reach a size of up to eight times the molecule. All mass spectra indicate only ions formed via Penning ionization induced by excited Argon molecules. As a result the most intense signal belongs to the $[M-H]^+$ molecule. To get further information about the oxidation process also the presence of water in the sample was analyzed. Depending on the amount of water the cluster size and the oxidation rate varies. Also the presence of atmospheric Oxygen and Nitrogen was investigated.

Conclusions

Measurements show that the main ionization product are the $[M-H]^+$ as well as $[M-2H-H]^+$. It is most likely that the corresponding mechanism is a radical reactions followed by a dissociation into a cation and one electron. This is also underlined by the formation of dehydrated species and C_xH_y losses. The cluster formation is highly dependent on structure of the analyte molecule.

Novel Aspect

Investigation of the clustering and fragmentation patterns of LTP ionization products under atmospheric conditions via FT-ICR-MS.

References

1. B. Raupers, Masterthesis, Christian-Albrechts-University Kiel, 2015.
2. M. Ihlenborg, B. Raupers, F. Gunzer, J. Grotemeyer, *Analyst*, 2015, 140, 7565. URL <http://dx.doi.org/10.1039/c5an01532d>.
3. A. W. Nørgaard, A. Vibenholt, M. Benassi, P. A. Clausen, P. Wolkoff, *Journal of The American Society for Mass Spectrometry* 2013, 24, 1090, URL <http://dx.doi.org/10.1007/s13361-013-0648-3>.

MP-239 / DEVELOPMENT AND APPLICATION OF A NOVEL PLASMA ION SOURCE WITH GC-MS/MS

Mehrnaz Sarrafzadeh - Reza Javahery - Charles Jolliffe

PerkinElmer, R&D, Woodbridge

Keywords: Tandem quadrupole mass spectrometry, Gas chromatography, Glow discharge plasma

Introduction

Developing an ion source for GC-MS/MS system capable of generating enhanced signal for molecular ions is very important for the analysis of unknown compounds. Traditional GC-MS with EI source often results in extensive fragmentation and the sample identification is mainly through a library search. Although library identification is a valuable tool, the reliability of the results cannot be ensured in many cases due to the absence of molecular ions.

Methods

All experiments were performed with a plasma ionization source integrated to a modified Qsight™ 220 triple quadrupole mass spectrometer equipped with a Clarus 680 GC. Several commercially available classes of compounds such as hydrocarbons, pesticides and terpenes were injected into the GC, which then were introduced to the plasma source by a heated transfer line, ionized and subsequently transmitted into the triple quad mass analyzer and detected.

Results

Our results indicated that the plasma source could efficiently ionize a wide range of compounds with a significant sensitivity enhancement of the precursor ions through soft ionization. The plasma source results in a very clean full mass spectrum with minimal fragments dominated by the precursor ions. This is in contrast to results of standard GC-EI systems in which the spectrum is dominated by fragments, while the precursor ion is absent or low in abundance. This allows for accurate identification of unknown compounds as well as direct quantitation by selection of precursor ions. Additionally, acquired fragmentation at various collision energies in MS/MS revealed very similar fragmentation patterns as those of electron ionization, which offers an opportunity for library identification.

Conclusions

Plasma source is a very simple, inexpensive and robust ion source which requires very low maintenance. GC-MS/MS with the plasma source offers improved performance and wider application over the conventional GC-MS or GC-MS/MS systems with the EI source by efficiently enhancing the sensitivity of GC-MS/MS analyses, with the unique advantage of producing high yield of intact molecular ions without extensive fragmentation due to the softer ionization.

Novel Aspect

Development of a low-cost and efficient plasma ionization source for a GC-MS/MS system which enhances molecular ions yield

MP-240 / A NOVEL VOLTAGE-ASSISTED PAPER SPRAY-APPI-HRMS METHOD FOR THE SCREENING AND QUANTIFICATION OF SEMI-VOLATILE FLUORINATED COMPOUNDS

Raquel Seró - Juan Francisco Ayala-Cabrera - Francisco Javier Santos - Encarnación Moyano

University of Barcelona, Chemical Engineering and Analytical Chemistry, Barcelona

Keywords: ambient mass spectrometry, paper spray, non-polar solvents, photoionization, semi-volatile fluorinated compounds

Introduction

Paper spray (PS) mass spectrometry is a fast ionization technique widely used for the detection of polar compounds mainly based in ESI ionization mechanisms. To improve the ionization of non-polar or low polar compounds, the traditional PS technique has been modified for the direct analysis of semi-volatile fluorinated pollutants by ambient-HRMS [1] in negative ionization mode using atmospheric pressure photoionization (APPI)[2].

Methods

Standards and samples were loaded into chromatography paper triangles and were placed in front of a Q-Orbitrap with the APPI-UV lamp set orthogonally between the sharp tip and the MS inlet. Toluene/cyclohexane was used for the extraction/transport of analytes to the paper tip and for dopant assisted APPI. By applying 1.5 kV, a spray was generated to desorb the non-polar analytes to be further photoionized in gas phase using negative ionization mode for HRMS analysis.

Results

Several working parameters, including dopant, solvent mixture, spray voltage and paper and APPI-UV lamp positions, which may affect the ionization, extraction and transport efficiency of the analytes, were carefully optimized for the analysis of the two most challenging families of fluorinated compounds by PS, zero-voltage PS-APPI and voltage-assisted PS-APPI. Under PS-HRMS and zero-voltage PS-APPI-HRMS conditions, most analytes were not or poorly ionized, while using voltage-assisted PS-APPI-HRMS, useful

$[M-H]^-$ and/or $[M+O_2]^-$ ions were detected. High resolution MS/MS experiments were also performed for further confirmatory purposes. Thus, major improvements in terms of sensitivity were obtained by using voltage-assisted PS-APPI ionization and method quality parameters demonstrated good repeatability, reproducibility and linearity with negligible matrix effects for the analysis of waterproof liquid matrices. Besides, internal standard method allowed the quantification of these families of compounds at pg mL^{-1} levels.

Conclusions

The use of APPI as an alternative ionization for PS-based techniques allowed the screening and quantification of non-polar and/or low polar compounds by PS. The voltage-assisted PS-APPI-HRMS method was applied for the detection of semi-volatile fluorinated compounds in real samples without any sample manipulation at $\mu\text{g mL}^{-1}$ levels demonstrating its potential as a fast and useful technique for quality control laboratories.

Novel Aspect

To the best of our knowledge, this is the first time that PS is successfully combined with APPI for the determination of semi-volatile fluorinated compounds.

References

1. Jahnke A., Berger U., J Chromatogr A, 1216, 410–421 (2009)
2. Kauppila T. J., Kostianen R., Anal Methods, 9, 4936–4953 (2017)

MP-241 / ADVANCES IN SELECTED ION FLOW-DRIFT TUBE MASS SPECTROMETRY, SIFDT-MS, FOR AMBIENT QUANTIFICATION OF TRACE VOLATILE COMPOUNDS IN AIR

Patrik Španěl - Anatolii Spesyvyi - David Smith - Kristýna Sovová

Czech Academy of Sciences, J. Heyrovský Institute of Physical Chemistry, Praha

Keywords: SIFT-MS, SIFDT-MS, VOCs, quantification, chemical ionization

Introduction

The Selected Ion Flow-Drift Tube Mass Spectrometry, SIFDT-MS, technique represents an important development of the Selected Ion Flow Tube Mass Spectrometry, SIFT-MS, technique building on its strengths (rapid switching between reagent ions and near thermal ion chemistry) and removing some of its weaknesses (lower sensitivity in comparison with proton transfer reaction, PTR, MS). It has a range of applications in biological, medical and environmental fields.

Methods

The selected reagent ions (either H_3O^+ , NO^+ or O_2^+) for selective chemical ionization of volatile compounds present in the air matrix are created in an external ion source and injected into the flow-drift tube reactor with well-defined reaction time and reagent ion-analyte molecule interaction energy, E_r . The absolute concentrations of the analytes can thus be accurately quantified using the chemical kinetics principles.

Results

Accuracy of quantification is shown to be determined by the accuracy of the reaction time measurement and by the uncertainties in the values of the ion-molecule reaction rate constants. The reactions of H_3O^+ , NO^+ and O_2^+ ions have thus been studied with a range of compounds, including ketones, terpenes, isoprene isoflurane, sevoflurane and acetonitrile in order to reveal their dependencies on ion-molecule interaction energies, E_r , whilst also gaining insight into their mechanisms. The rate constants, k , for the reactions of

H_3O^+ proceeding via exothermic proton transfer are seen to reduce by about one third over the Er range from 0.05 eV up to 0.5 eV, in agreement with the collisional theory. The NO^+ /ketones reactions often proceed via the formation of excited $(\text{NO}+\text{M})^*$ adduct ions that partially fragment, and k reduce with Er some times more dramatically. The k for the O_2^{++} charge transfer reactions are observed to be greater than their calculated kc implying that long distance electron transfer occurs.

Conclusions

Whilst several advantages of SIFDT-MS over the established SIFT-MS are becoming apparent (including better sensitivity, improved accuracy and feasibility of design of portable instrument), it is very important to assess the potential complications due to increased Er, including fragmentation and changes of k. This technique thus has a potential for accurate analyses in matrices such as environmental air, food flavours or exhaled breath.

Novel Aspect

New experimental data are presented and interpreted on the kinetics of ion-molecule reactions allowing a new approach to soft chemical ionization of trace compounds in air.

References

1. Spesyvyi A., Smith D., Španěl P., *Analytical Chemistry*, 87, 12151 (2015)
2. Spesyvyi A., Sovová K., Španěl P., *Rapid Communications in Mass Spectrometry*, 30, 2009 (2016)
3. Spesyvyi A., Smith D., Španěl P. *Phys. Chem. Chem. Phys.* 19, 31714(2017)

MP-242 / A TRIPLE REAGENT ION SOURCE FOR TRACE VOCs QUANTIFICATION BY SELECTED ION FLOW-DRIFT TUBE MASS SPECTROMETRY

Anatolii Spesyvyi - Patrik Španěl

J. Heyrovský Institute of Physical Chemistry of the CAS, Department of Chemistry of Ions in Gaseous Phase, Prague

Keywords: SCI-MS, SIFDT-MS, VOCs, quantification, chemical ionization

Introduction

Selected Ion Flow-Drift Tube Mass Spectrometry, SIFDT-MS, was developed to quantify volatile organic compounds in air. The original instrument had two quadrupole mass filters: to select the reagent ions from the microwave discharge into the drift tube reactor and to analyze the product ions. To increase the sensitivity, the first quadrupole and the microwave discharge were replaced by a triple source creating the reagent ions by sequences of ion-molecule reactions.

Methods

The SIFDT-MS apparatus and the basic principles of the technique were previously described [1, 2]. For the present work a newly designed ion source was constructed consisting of three hollow cathode discharges arranged parallel to the drift tube axis. Different gas (water vapor, NO or O_2) is used in each discharge creating reagent ions H_3O^+ , NO^+ or O_2^+ rapidly switchable by drift lens voltages.

Results

It was demonstrated that the ion source allows simultaneous use of three different gases to produce reagent ions of interest that can be selected by applying optimized sets of DC electric potentials to drift lenses. Increase in the analytical sensitivity is described that was achieved by greater reagent ion intensity increased several times due to the absence of the orifice between the ion source and the drift tube, which are now pumped by a common vacuum line and have similar pressures of about 2 mbar. Ion-molecule

reactions in the drift tube are governed by the voltage applied across the drift tube. It is shown how the relative electric field strength E/N affects the reaction time, the ion interaction energy and the diffusion losses. Illustrative results include enhanced product ion fragmentation [3] and suppression of unwanted reactions, for example formation of hydrated ions.

Conclusions

The presented enhancement of SIFDT-MS makes this technique more sensitive and allows construction of robust and small instruments. It thus represents a viable addition to the family of soft chemical ionization mass spectrometric techniques such as SIFT-MS and PTR-MS that are becoming widely used in the area of real time quantification of the VOCs in the medical, food and environmental fields.

Novel Aspect

Three reagent ions H_3O^+ , NO^+ or O_2^+ are produced in the triple ion source simultaneously and can be selected instantly at greater intensities than possible using a mass filter.

References

1. Spesyvyi A., Smith D., Španěl P., *Analytical Chemistry*, 87, 12151 (2015)
2. Spesyvyi A., Španěl P., *Rapid Communications in Mass Spectrometry*, 29, 1563 (2015)
3. Spesyvyi A., Sovová K., Španěl P., *Rapid Communications in Mass Spectrometry*, 30, 2009 (2016)

MP-243 / AIRBORNE LASER-SPARK ION SOURCE FOR DIRECT MICROFLUIDIC COUPLING

Sebastian van Wasen⁽¹⁾ - ***Andreas Bierstedt***⁽²⁾ - ***Jens Riedel***⁽²⁾ - ***Dietrich A. Volmer***⁽¹⁾

Humboldt-Universität zu Berlin, Department of Chemistry, Berlin (1) - Federal Institute for Materials Research and Testing, BAM, Berlin (2)

Keywords: ambient mass spectrometry, laser-induced plasma, DPSS laser, high repetition rate

Introduction

The development and improvement of new ionization techniques for mass spectrometry often requires dedicated, specific sampling approaches. Recently, a novel ionization scheme for ambient MS has been introduced based on a quasi-continuous laser-induced plasma, which was ignited directly before the MS inlet. This setup combines the general advantages of ambient ionization, provides electro neutrality, sufficient duty cycle and a ubiquitous plasma medium. [1, 2]

Methods

A high repetition rate DPSS laser (Conqueror 3-LAMBDA, Nd:YVO4, 1 - 500 kHz, average output power: 12 W at 50 kHz, Compact Laser Solutions GmbH, Germany) and the corresponding optomechanical system were installed on an optical breadboard above the inlet of a LCQ DecaXP ion trap mass spectrometer. The quasi-continuous airborne plasma was ignited inside the sprayed sample in front of the inlet via focused laser irradiation.

Results

The introduction of liquid samples into laser-induced plasmas requires higher plasma power during solvent evaporation as compared to gaseous samples. This increased demand was approached via a two-fold strategy: Firstly, an alternative, more powerful, laser plasma driven by the fundamental instead of the second harmonic wavelength was implemented, which provided a 10-fold increase of signal intensity, while maintaining the same reagent ion pattern as the previous plasma. Protonated water clusters $[(H_2O)_nH]^+$, NH_4^+ as well as charge transfer promoting ion O_2^+ , dominated the reagent ion mass spectrum. Secondly, a miniaturized nebulizer was used to minimize the size of the plasma quenching solvent droplets. The result

of these improvements was a new and very stable ion source for direct microfluidic coupling. A variety of samples demonstrated the performance of the ion source.

Conclusions

A laser-driven plasma was shown to be a powerful ion source for gaseous and solid samples. For the first time, liquid samples were examined using the novel source. In addition to demonstrating an improved strategy for igniting the laser plasma, this contribution also covers the miniaturization of the spray source for enhanced ionization, while minimizing sample consumption via a microfluidic spray systems.

Novel Aspect

The first laser plasma ionization spectrum was recorded for liquid sample introduction.

References

1. Bierstedt, A. ; Riedel, J., *Eur. J. Mass Spectrom.* 2016, 22(3), 105-114 (2016)
2. Bierstedt, A. et al., *Anal. Chem.*, 89(6), 3437-3444 (2017)

MP-244 / DEVELOPMENT AND THE APPLICATION OF A 96 HIGH THROUGHPUT SOLID PHASE MICROEXTRACTION (SPME) TRANSMISSION MODE (TM) SYSTEM FOR DIRECT ANALYSIS IN REAL TIME (DART) MASS SPECTROMETRY (MS)

Tijana Vasiljevic - Germán Augusto Gómez-Ríos - Janusz Pawliszyn

University of Waterloo, 200 University Avenue West, Waterloo

Keywords: High throughput; Ambient mass spectrometry (AMS); Direct analysis in real time (DART); Solid phase microextraction (SPME); Drugs

Introduction

Development of a high throughput system capable of delivering the two “Rs” (rapidity and reliability) is one of the most sought-after research fields in the area of ambient ionization mass spectrometry (AIMS). In this research, we have taken the single SPME-TM mesh and transformed it into a 12 mesh format for rapid high throughput detection/quantification of opioids in biomatrices using direct analysis in real time mass spectrometry (MS).

Methods

SPME-TM meshes were manufactured to a 12 mesh format using pre-made photochemically milled stainless steel substrates, which were etched and adequately coated with C18 particles. Urine and plasma samples were spiked with a set of opioids, followed by a 10 minute extraction step and desorption at 450°C of DART source with tandem MS analysis. A set of 96 samples was extracted, analyzed and interpreted in less than 1.5 hours.

Results

The SPME-TM 12 mesh strip were created utilizing a careful fabrication protocol with respect to etching and coating process, while their performance was maximized using a robust holder made at University of Waterloo (UW) capable of analyzing all 12 meshes at the same time. Positive preliminary assessment of mesh performance (*i.e.* coating stability, extraction efficiency) led us to apply the meshes for a high throughput quantification of selected opioids in urine and plasma. Experimental procedure was optimized with respect to extraction time, DART source temperature (TEM), rail speed and pH of samples ultimately resulting in low limits of quantification (LOQ) (*i.e.* 0.5 ng mL⁻¹ fentanyl for PBS and urine and 1 ng mL⁻¹ for plasma), acceptable accuracy (73.4-109.5 %) and repeatability (1.6-11.3%). Indeed, these results take a firm

step toward establishing DART SPME-TM as potential tool for detection of varying substances in scenarios where high throughput is required.

Novel Aspect

First reported use of a high throughput SPME 96 system for opioid detection using DART/MS

MP-245 / SINGLE-USE POLY(ETHERETHERKETONE) SOLID-PHASE MICROEXTRACTION-TRANSMISSION MODE DEVICES FOR RAPID SCREENING AND QUANTITATION OF DRUGS OF ABUSE IN ORAL FLUID AND URINE VIA DIRECT ANALYSIS IN REAL-TIME TANDEM MASS SPECTROMETRY

Tijana Vasiljevic - Germán Augusto Gómez-Ríos - Janusz Pawliszyn

University of Waterloo, Chemistry, Waterloo

Keywords: Direct-analysis in real time (DART), Solid phase microextraction (SPME), Polyetheretherketone (PEEK), Drugs, Quantification

Introduction

Rapid detection of prohibited substances in biologicals is important in forensic science. Herein, we report on results obtained for detection and quantification of drugs of abuse in urine and oral fluid (OF) by using biocompatible single use plastic SPME-TM meshes made of polyetheretherketone (PEEK) coupled to direct analysis in real time mass spectrometry (DART-MS).

Methods

PEEK meshes were coated with hydrophilic lipophilic (HLB) particles. The biofluids were spiked with a mix of drugs and internal standards. The meshes were used for extraction of selected drugs in 700 μ L of urine and 300 μ L of OF. Extraction time was kept at 1 minute, followed by wash in water and drying. Meshes were placed onto a homemade holder and desorbed between the DART source and the MS on a linear rail. Analysis was carried out using a triple quadrupole MS with MRM mode.

Results

The goal behind development of the plastic meshes was to be create a tool which could potentially be used for on-site in vivo analysis. Indeed, the plastic meshes provided us good results for quantification of 13 different drugs in urine and OF, with linearity ranging from 0.5-200 ng/mL and R2 values > 0.99. Validation points were tested at 4 different levels and provided acceptable accuracy in the range of 75- 127 %. Also, the LOQs obtained for the analytes were below the limits imposed by substance abuse and mental health services administration (SAMSHA). In addition, OF samples were collected from a volunteer for detection of caffeine levels prior to and after coffee consumption. The results obtained at different time points were in accordance with caffeine half-life for healthy adults. In addition, we looked at impact of material on heat transfer in DART. Such findings highlight the fact that ambient ionization techniques such as DART are very useful tools in providing law enforcement with much needed reliable and rapid data.

Conclusions

Results obtained indicate great potential of SPME-TM DART for use in many scenarios where rapid, reliable and quality data is required in a timely fashion. Utilization of alternate materials such as PEEK opens up a great prospect for use in situations where in vivo detection of prohibited substances may be required, in addition to having great prospect for use with portable MS instrumentation.

Novel Aspect

First use of biocompatible Poly(etheretherketone) as SPME-TM meshes for rapid detection of drugs of abuse and exploration of impact of material on heat transfer in DART

References

1. Vasiljevic, T. ; Gómez-Ríos, G. A. ; Pawliszyn, J. Analytical Chemistry 2017, 90 (1), 952

MP-246 / IMPROVEMENTS IN SMALL MOLECULE DETECTION VIA MALDI-MS USING BINARY MIXTURES OF NOVEL IONLESS MATRICES

Jerrit Weissflog - Aleš Svatoš

Max-Planck-Institute for Chemical Ecology, Research Group Mass Spectrometry/Proteomics, Jena

Keywords: MALDI-MS, ionless matrices, binary matrices, small molecule detection

Introduction

Matrix-assisted laser desorption/ionisation mass spectrometry (MALDI-MS) is mainly used for the analyses of larger molecules ($m/z > 1000$) due to interfering mass signals below m/z 600 generated by conventional matrices. So called ionless matrices like graphene [1] and the recently introduced DPN [2] do not generate matrix signals and therefore allow detection of small molecules. This study aims to improve their sensitivity by applying mixtures of such matrices.

Methods

Several established matrix compounds like 1,8-di(piperdiny)-naphthalene[DPN], 9-aminoacridine[3] and 2-naphthylsulfonic acid[4] were mixed in distinct ratios and their ionization potential was assessed via MALDI-ToF-MS using a variety of analytes *e.g.* fatty acids and amino acids. Additionally, in several experiments graphene, as well as strong acids or bases were applied together with these matrices to test their influence on the ionization performance.

Results

When the ionless matrices are used in binary mixtures, a clear increase in signal intensity can be observed. The ionizing properties of these blends are higher than those of the single matrices. Graphene for example was able to improve the performance of DPN in negative ion mode especially for hydrophobic analytes like fatty acids while the spectra remained free of any signals besides the ones from the analytes. Furthermore, the required energy per laser pulse for optimal signal intensity is reduced when graphene is applied along side DPN.

Conclusions

The mechanism of ionization in MALDI is a complex interplay between absorption and distribution of laser energy, proton transfer between analyte and matrix, and separation of ion pairs in the gas phase. By employing mixtures of matrices, a possible existing shortcoming of a single matrix in one of these aspects can be mitigated if the other matrix surpasses the first in this particular property.

Novel Aspect

A substantial improvement of sensitivity for the detection of small molecules via MALDI-MS using mixtures of ionless MALDI-matrices as well as graphene.

References

1. Dong X., Cheng J., Li J., Wang Y., Anal. Chem., 82, 6208 (2018).
2. Weißflog J., Svatoš A., RSC Adv., 6, 75073 (2016).
3. Miura D., Fujimura Y., Yamato M., Hyodo H., Tachibana H., Wariishi H., Anal. Chem., 82, 9789 (2010).
4. Shroff R., Rulisek L., Doubsky J. Svatoš A., PNAS, 106, 10092 (2009)

MP-247 / MATRIX-ASSISTED-IONIZATION-IN-VACUUM COUPLED TO HIGH-RESOLUTION FOURIER TRANSFORM MASS SPECTROMETRY (MAI-FTMS) IN CHARACTERIZING LIPIDS IN COOKING OIL AND BACTERIA

Rohanna Liyanage⁽¹⁾ - *Jennifer Gidden*⁽¹⁾ - *Jackson Lay, Jr.*⁽²⁾ - *Charles Wilkins*⁽²⁾

University of Arkansas, Statewide Mass Spectrometry Facility, Fayetteville, United States (1) - University of Arkansas, Chemistry and Biochemistry, Fayetteville, United States (2)

Keywords: MAI, Ionization, Matrix, Lipids, Bacteria

Introduction

MAI seems to be a reliable alternative to ESI and MALDI for some applications. Solvent free and laser free aspects of MAI appeal to most mass spectrometrists [1,2]. The primary disadvantage of MAI is rapid evaporation/consumption of the matrix in the vacuum during analysis. There are relatively few literature reports of measuring lipids by MAI-FTMS. Here we show the use of MAI-FTMS lipid profiles in differentiating a variety of cooking oils and bacteria.

Methods

Samples were prepared by dissolving in matrix solutions. MAI analysis was done by FTMS immediately after drying the sample solutions in either positive or negative ion mode. The FT-ICR mass spectrometer is a 9.4 T Apex FT-ICR mass spectrometer (Bruker, Bremen, Germany). After air drying and loading, the plate voltage was set at 300 V. The acquisition size was at least 512 K points for higher signal abundance. Twenty acquisitions were averaged to obtain a mass spectrum.

Results

3-nitrobenzonitrile (3-NBN) solubilized in acetonitrile was used as matrix. Quality MAI-FTMS spectra from TAGs in cooking oil result when prepared with high concentrations of matrix dissolved in pure acetonitrile. Phospholipids from the bacteria only required dilute matrix. Generally, increasing the amount of water in the matrix allowed high quality phospholipid spectra with increases in the lifetime of matrix in vacuum. MS/MS fragmentation yields all important peaks in MAI-FTMS. TAGs profiles obtained for all the cooking oils, peanut, soy, sunflower, coconut, cocoa butter, canola, and vegetable, are identical to those obtained with MALDI. MAI had more reproducible intensity variation, allowing statistical determination of small compositional differences [3]. Although most of the MALDI-observed phosphatidylethanolamines (PE) lipids from *E. coli* are adducts of sodium or potassium, MAI primarily produces protonated

Conclusions

The newly-developed MAI-ionization technique is a rapid and simple method for characterization of polar and non-polar lipids in complex mixtures [1,2]. This technique is complementary to ESI and MALDI mass spectrometry methods and combines some of the best features of both.

Novel Aspect

Application of MAI-FTMS in characterizing polar and non-polar lipids in complex samples

References

1. Trimpin S. J. Am. Soc. Mass Spectrom. "Magic" Ionization Mass Spectrometry 2016; 27: 4-21.
2. Wang B, Tisdale E, Trimpin S, Wilkins CL. Matrix-Assisted Ionization Vacuum for High-Resolution Fourier Transform Ion Cyclotron Resonance Mass Spectrometers. Anal. Chem. 2014; 86: 6792-6796.
3. Lay JO, Gidden J, Liyanage R, Emerson B, Durham B. Rapid characterization of lipids by MALDI MS. Part 1: Bacterial taxonomy and analysis of food oils. Lipid Technology. Lipid Technology: 24 (1) 11-14
4. Lay JO, Liyanage R, Durham B, Brooks J. Rapid characterization of edible oils by direct matrix-assisted laser desorption/ionization time-of-flight mass spectrometry analysis using triacylglycerol's. Rapid Commun. Mass Spectrom. 2006, 20, 952-958.

MP-248 / 3D-PRINTED LOW-TEMPERATURE PLASMA (3D-LTP) PROBE FOR IMAGING AND IN VIVO STUDIES

Sandra Martínez-Jarquín⁽¹⁾ - *Abigail Moreno-Pedraza*⁽¹⁾ - *Héctor Guillén-Alonso*⁽¹⁾ - *Humberto Herrera-Ubaldo*⁽²⁾ - *Ignacio Rosas-Román*⁽¹⁾ - *Cesaré Ovando-Vázquez*⁽³⁾ - *Stefan de Folter*⁽²⁾ - *Robert Winkler*⁽¹⁾

CINVESTAV Unidad Irapuato, Biotechnology and Biochemistry, Irapuato (1) - CINVESTAV UGA, LANGEBIO, Irapuato (2) - IPICYT, CNS, San Luis Potosi (3)

Keywords: low-temperature plasma, 3D-printing, ambient imaging, in vivo studies, real-time monitoring

Introduction

A wide range of organic molecules from solid, liquid and gaseous samples can be detected with low-temperature plasma (LTP) ionization mass spectrometry (MS) under ambient conditions [1]. LTP probes are relatively easy to build, but there is no commercial system available. Thus, we designed an LTP probe from 3D-printed and standard parts and tested its analytical performance on chemical and biological samples.

Methods

We constructed an LTP probe from standard and 3D-printed components and tested the suitability of the 3D-LTP probe for MS by measuring reference standards and complex samples[2]. For automated sampling and imaging, we built a robot based on RepRap technology. We investigated the biocompatibility of LTP-MS by monitoring the nicotine content in tobacco leaves for 10 days and subsequently evaluating the tissue damages with a microscope.

Results

The 3D-LTP probe was functional and coupled to a standard MS analyzer. The 3D files were published under a Creative Commons (non-commercial) Public License.

Compounds of different volatility and polarity and with a molecular weight up to ~500 Da could be detected from simple and complex matrices. The diameter of the plasma beam was about 200 μm . The technical resolution of the sampling robot was 12.5 μm . Therefore, ambient mass imaging of eukaryotic tissues with cellular resolution is possible when combined with laser desorption (LD).

We subjected tobacco plants with modulated nicotine biosynthesis to daily measurements for 11 days. The LTP-MS measurements resembled the expected time course which indicated correct monitoring of the biosynthesis. Staining and microscopic evaluation only showed minor tissue damage and no systemic effect was observed. This evidence demonstrated the high potential of LTP-MS for in vivo biosynthesis studies[3].

Conclusions

The 3D-LTP probe can be easily built from 3D-printed and commercially available parts. The plasma jet diameter is sufficient for imaging macroscopic samples. Short-time exposure of plant surfaces to LTP only resulted in negligible (microscopic) lesions, thus demonstrating the value of LTP-MS for in vivo measurements. Our license model promotes the free exchange and innovation of LTP probe designs.

Novel Aspect

The Creative Commons (non-commercial) license of the 3D-LTP probe permits the rapid building of custom low-temperature plasma probes.

References

1. Martínez-Jarquín, S. & Winkler, R., *TrAC Trends Anal. Chem.*, 89, 133–145 (2017).
2. Martínez-Jarquín, S., Moreno-Pedraza, A., Guillén-Alonso, H., & Winkler, R., *Anal. Chem.*, 88, 6976–6980 (2016).
3. Martínez-Jarquín, S., Herrera-Ubaldo, H., de Folter, S., & Winkler, R., *Talanta*, 185, 324–327 (2018).

MP-249 / IDENTIFICATION OF MEDULLOBLASTOMA SUBGROUPS ON INTRASURGICAL TIMESCALES WITH PICOSECOND INFRARED LASER DESORPTION MASS SPECTROMETRY

Michael Woolman⁽¹⁾ - *Claudia Kuzan-Fischer*⁽²⁾ - *Isabelle Ferry*⁽²⁾ - *Taira Kiyota*⁽³⁾ - *Betty Luu*⁽²⁾ - *Megan Wu*⁽²⁾ - *David Munoz*⁽⁴⁾ - *Sunit Das*⁽²⁾ - *Ahmed Aman*⁽³⁾ - *Michael Taylor*⁽²⁾ - *James Rutka*⁽²⁾ - *Howard Ginsberg*⁽⁴⁾ - *Arash Zarrine-Afsar*⁽⁵⁾

University of Toronto, Department of Medical Biophysics, Toronto (1) - Hospital for Sick Children, University of Toronto, Toronto (2) - Ontario Institute for Cancer Research, University of Toronto, Toronto (3) - St Michael's Hospital, University of Toronto, Toronto (4) - University Health Network, University of Toronto, Toronto (5)

Keywords: Medulloblastoma, picosecond infrared laser, mass spectrometry, molecular pathology, real time analysis

Introduction

Medulloblastoma (MB) is the most common brain tumor in children that is comprised of four different subgroups (WNT, SHH, Group 3, Group 4) each with distinct prognostic characteristics [1-3]. While risk stratification of MB patients based on molecular features may offer the prospect of personalized therapies, the long analysis times of conventional classification methods prevents subgroup-specific adjustments of the extent or the aggressiveness of the resection.

Methods

We investigated the potential of rapid tumor small molecule profiling with Picosecond Infrared Laser desorption Mass Spectrometry (PIRL-MS) [4,5] as a new method of medulloblastoma subgroup classification. A local tissue bank comprising 113 ex vivo medulloblastoma tumors was subjected to 10-second PIRL-MS data collection, and Principal Component Analysis with Linear Discriminant Analysis (PCA-LDA).

Results

A model for MB subgroups was built based on the laser extracted lipid profiles of 72 independent tumors, comprised of 190, 10-second PIRL-MS data points. Forty one de-identified unknown tumors were subjected to 10-second PIRL-MS sampling and real time PCA-LDA analysis using the model above where the resultant 124 sampling events, after the application of a 95% PCA-LDA prediction probability threshold, yielded a 99% correct classification rate, and only one misclassified data point. This was accompanied with 28 sampling events that were not classifiable into any known MB subgroups. We thus used an extended MB signature library that was comprised of 278 sampling events to better capture the intrinsic biological heterogeneity of MB, through which we were able to correctly reclassify 21 out of the said previously unclassifiable sampling events with no misclassification. Detailed post ablation histopathological analysis of MB samples suggests that intratumoral heterogeneity at the site of laser ablation can explain the apparent misclassification.

Conclusions

Our observations suggest that rapid tumor small molecule content profiling with PIRL-MS is effective in medulloblastoma subgroup classification with only 10 seconds of sampling, an analysis time thus opening up the possibility of subgroup specific treatment of medulloblastoma beginning in the operating room during a surgical procedure.

Novel Aspect

Rapid intraoperative identification of Medulloblastoma subgroups in 10 seconds based on lipid profiling with Picosecond InfraRed Laser Mass Spectrometry is now possible.

References

1. Northcott, P. A. et al. Medulloblastoma comprises four distinct molecular variants. *J Clin Oncol* 29, 1408-1414 (2011).
2. Ostrom, Q. T. et al. CBTRUS Statistical Report: Primary Brain and Central Nervous System Tumors Diagnosed in the United States in 2008-2012. *Neuro Oncol* 17 Suppl 4, iv1-iv62, (2015).
3. Ostrom, Q. T. et al. Alex's Lemonade Stand Foundation Infant and Childhood Primary Brain and Central Nervous System Tumors Diagnosed in the United States in 2007-2011. *Neuro Oncol* 16 Suppl 10, x1-x36 (2015).
4. Woolman, M. et al. Optimized Mass Spectrometry Analysis Workflow with Polarimetric Guidance for ex vivo and in situ Sampling of Biological Tissues. *Scientific reports* 7(1): 468-480 (2017).
5. Woolman, M. et al. Rapid determination of medulloblastoma subgroup affiliation with mass spectrometry using a handheld picosecond infrared laser desorption probe. *Chemical Science* 8, 6508-6519 (2017).

MP-250 / IMPROVEMENT OF NANO ELECTROSPRAY AND DESORPTION NANO ELECTROSPRAY SIGNAL BY REGULATION OF PRESSURE IN AN INLET OF TRIPLE QUADRUPOLE MASS SPECTROMETER

Lucie Borovcová⁽¹⁾ - Petr Fryčák⁽¹⁾ - Dominik Mosler⁽¹⁾ - Štěpán Dostál⁽¹⁾ - Vladimír Havlíček⁽²⁾ - Karel Lemr⁽¹⁾

Palacký University, Department of Analytical Chemistry, Olomouc, Czech Republic (1) - ASCR, Institute of Microbiology, Prague, Czech Republic (2)

Keywords: mass spectrometry, nanoelectrospray ionization, desorption nanoelectrospray

Introduction

Both nanoelectrospray (nanoESI) and desorption nanoelectrospray (nanoDESI) [1,2] use a narrow spray tip and consume small amount of samples. We modified a commercial ion source to switch between both ionization modes easily by the change of direction of a spray but sensitivity of nanoDESI was not sufficient. To improve ion signal intensity we have tested regulation of pressure drop between an atmospheric and evacuated region of a mass spectrometer.

Methods

Experiments were carried out using a Xevo TQD triple quadrupole mass spectrometer (Waters, Manchester, UK). The commercial inlet was replaced by a custom-made inlet with a sample cone and a heated capillary allowing for the pressure regulation. Nanoelectrospray tips ($2 \pm 1 \mu\text{m}$ I. D., PicoTips, New Objective, Woburn, USA) were used for nanoESI as well as nanoDESI. Compounds in the wide mass range (m/z from 111 to 1922) were ionized to test the performance of ion sources.

Results

Mass spectra were acquired for different values of pressure drop between an atmospheric and evacuated region and absolute as well as relative ion signal intensities were compared. Repeatability of the intensities was also evaluated. The relative intensities of ions of analytes were not influenced significantly but absolute signal was changed by pressure regulation in the whole studied mass range. Using nanoESI in positive ion mode, absolute ion signal of oxidized ultramark 443 (m/z 459) and ultramark 1621 (m/z 1222) increased 5 times for pressure drop 0.4 bar. For nanoDESI and the same pressure drop, signal of ion at m/z 459 and m/z 1222 was enhanced 7 and 9 times, respectively. Similar behavior was also observed in negative ion mode, e.g. for benzoic acid, sorbic acid or acetylsalicylic acid. Relative standard deviations of signal intensities did not change significantly for different values of pressure drop.

Conclusions

New design of the mass spectrometer inlet allowed to increase yield of ions. For nanoDESI some ions were

not detected without the pressure regulation. The pressure drop between 0.4-0.7 bar was suitable for the whole studied m/z range. The source modification needed an additional rough pump but sensitivity was improved dramatically. The effect can be influenced by a pumping system of a mass spectrometer as the pressure in an evacuated region depends on its capacity.

Novel Aspect

New design of an ion source of a triple quadrupole mass spectrometer allowing pressure drop regulation to increase ion signal intensity.

References

1. V. Ranc, V. Havlíček, P. Bednář, K. Lemr; Chem. Listy, 101, 524 (2007).
2. L. Hartmanová, I. Lorencová, M. Volný, P. Fryčák, V. Havlíček, H. Chmelíčková, T. Ingr, K. Lemr; Analyst, 141, 2150 - 2154 (2016).

Acknowledgment

The authors gratefully acknowledge the support from the Czech Science Foundation (Grant 16-20229S).

MP-251 / ACOUSTIC MIST IONISATION PLATFORM FOR HIGH-THROUGHPUT MS-BASED DRUG DISCOVERY

Phillip DeLand⁽¹⁾ - *Ian Sinclair*⁽²⁾ - *Gareth Davies*⁽³⁾ - *David Murray*⁽³⁾ - *Elizabeth Mouchet*⁽³⁾ - *Martin Bachman*⁽⁴⁾ - *Mattias Rohman*⁽⁵⁾ - *Jon Wingfield*⁽⁶⁾

Labcyte Inc., Mass Spectrometry Product Management, San Jose, United States (1) - AstraZeneca, Sample Management, Alderley Park, United Kingdom (2) - AstraZeneca, Hit Discovery, Alderley Park, United Kingdom (3) - Medicines Discovery Catapult, Drug Discovery, Alderley Park, United Kingdom (4) - AstraZeneca, Discovery Biology, Molndal, Sweden (5) - AstraZeneca, Mechanistic Biology & Profiling, Cambridge, United Kingdom (6)

Keywords: acoustic loading; mass spectrometry; high-throughput screening;

Introduction

Mass spectrometry (MS) is a widely utilized label-free technology capable of measuring a broad range of analytes but lacks the throughput to support large HTS campaigns. Current commercial ESI MS-based systems using solid-phase extraction (SPE) are capable of running medium-sized screens. We have implemented a platform utilizing acoustics to load samples into a time-of-flight mass spectrometer capable of analysing more than 100,000 samples per day.

Methods

An acoustic transducer emits sound waves into a sample in a 384-well plate located on a moving XY-stage. High voltage is applied to a charging cone above the well, inducing charge separation in the sample. The acoustic wave ejects a mist of droplets which are drawn through a heated transfer interface into the MS, during transit droplets dry and an ion beam is generated. The extracted ion peak for substrate and product from each well are used to calculate enzyme activity.

Results

AMI-MS was used to build an enzyme assay to identify inhibitors from the AZ compound collection. Having established the kinetic parameters for this enzyme target, validation sets were tested to understand the variability within the assay. Once the robust nature of both the assay and this new MS technology were proven to meet the requirements of a high throughput screen the assay was moved into production. In total over 2 million samples were tested using a single automated AMI-MS platform over a 5 week period. Typical throughput of 100,000 samples per day were achieved. From the collection, 6745 compounds were identified in the first round of screening as actives - a hit rate of 0.3%. From this population of compounds, 4872 were followed up in concentration response format to confirm the activity, ~75% were confirmed as

actives. The AMI-MS platform was further utilised during triaging of compounds to identify false positive samples which were likely to be contaminated with zinc.

Conclusions

AMI-MS proved a robust technology capable of screening a large compound collection (>2 million samples) and identifying chemical equity capable of inhibiting the enzyme target. Typically, this single instrument was capable of processing samples at a rate of 100,000 data points per day. The final hit rate in this assay was ~0.3%, of these the majority confirmed as active in follow up testing (75%).

Novel Aspect

This is the first demonstration of AMI-MS technology successfully supporting a full collection screen against a pharmaceutically relevant target.

MP-252 / STUDY ON THE IONIZATION MECHANISM OF POLYMERS UNDER CLUSTER ION IRRADIATION

Makiko Fujii⁽¹⁾ - ***Jiro Matsuo***⁽²⁾

Yokohama National University, Graduate School of Environment and Information Sciences, Yokohama, Japan (1) - Kyoto University, Quantum Science and Engineering Center, Kyoto, Japan (2)

Keywords: ionization, cluster ion beam, polymer, SIMS

Introduction

Secondary Ion Mass Spectrometry (SIMS) has been widely used for the analysis of semiconductors and electronic devices in past decades. Recently, the development of novel ion beam, cluster ion beam, enabled us to detect large organic molecule in SIMS analysis. However, the sensitivity of molecular ion derived from large organic molecules are still insufficient for high-resolution imaging mass spectrometry with SIMS.

Methods

To increase the sensitivity of organic molecules in SIMS, some alkali metals and organic acids were added to the sample. This is because molecular ions are generally detected as alkali metal-adducted molecular ion or proton-adducted molecular ions. Therefore, by adding alkali metal ion or proton, the molecular ion yields can be thought to be enhanced.

Results

The model polymer samples used in this study were Polyethylene Glycol (PEG) of average molecular weight of 1000, 2000 and 3000. Lithium, Sodium, Potassium and Cesium were added to PEG sample to investigate the molecular ion yield enhancement by cationization agent. Citric acid and some other organic acids were added to investigate the molecular ion yield enhancement by protonation agent. The concentration of each additive was optimized in the preliminary study. The ion beam used in this study was Argon gas cluster ion beam (Ar-GCIB) of acceleration energy of 10 keV originally developed in Kyoto University. The ion dose density was around $1 \text{ E}+12 \text{ ions/cm}^2$ in each measurement. The molecular ion yield of PEG1000, 2000 and 3000 were calculated and compared each other in order to discuss the degree of the molecular ion yield enhancement depends on the molecular weight.

Conclusions

The molecular ion yield enhancement could be confirmed in all additives in the comparison with without additive sample. In addition, the degree of the molecular ion yield enhancement depended on the chemical structure of targeted organic molecule besides molecular weight.

Novel Aspect

The molecular ion yield with cationization or protonation agents was investigated. It will help

understanding the ionization mechanism under ion beam irradiation.

MP-253 / A UNIQUE APPROACH TO CHARGING IONS FOR LCMS

Charles Jolliffe - Lisa Cousins - Reza Javahery - Anna Kornilova - Heather Gamble - Mehrnaz Sarrafzadeh

PerkinElmer Inc., RND, Woodbridge

Keywords: ESI, APCI, MS/MS, LCMS

Introduction

This paper discusses a novel approach to charging ions at atmospheric pressure, suitable for both polar and nonpolar compounds, as well as extracting ions from opposite polarity droplets. Typically ESI is used with the majority of polar compounds and APCI with less polar. Here we discuss a single approach capable of producing an optimum signal for the full range of analyte polarity and to extract ions of opposite polarity.

Methods

Experiments were performed using a QSight 210 triple quadrupole mass spectrometer with a modified ion source probe that permits adjustable regions of turbulence and laminarization at the exit of the ion formation region. Analytes including myoglobin, reserpine and chlordane were studied. MRM sensitivity was taken at a 500 μ l/m, as a function of droplet charging. The same single probe physical position and sample tube position was maintained.

Results

Results are reported for a prototype ion source probe that embodies the characteristics of a turbulence and ion entrainment as well as the ability to optimally charge droplets, including the selection of an ion polarity opposite to the droplet polarity. Single probe S/N at high flow is compared with typical individual APCI and ESI probe results, using compounds which respond to both ESI and APCI, as well to those which respond predominantly to ESI or APCI. For example, although conventionally myoglobin is ESI favourable, reserpine is responsive to both ESI and APCI, and chlordane is APCI favourable, the new probe demonstrates suitable performance for all types. Single probe MRM data was taken of the opposite polarity of the charged droplets for all of the chosen analytes. For example, it was observed that the reserpine droplets can have a variable negative charge while still observing positive reserpine ions. A graph of the reserpine sensitivity vs voltage difference between the probe nebulizer tip and desolvation region illustrates this phenomenon.

Conclusions

The data from the MRM of the chosen analytes show that single probe sensitivity with optimized droplet charging is at least comparable to the individual APCI or ESI probes.

In addition, the single probe shows that significant MRM sensitivity is possible for APCI favourable compound chlordane using oppositely charged droplets. The effect of oppositely charged droplets upon chemical interferences is unexplored at this juncture.

Novel Aspect

A single atmospheric pressure ion probe produces comparable sensitivity to individual ESI and APCI probes as well as producing ions from oppositely charged droplets.

MP-254 / DEVELOPMENT OF MASS SPECTROMETRY-BASED EXPLOSIVES TRACE DETECTORS FOR SECURE CITY

Shun Kumano⁽¹⁾ - *Masuyuki Sugiyama*⁽¹⁾ - *Yasuaki Takada*⁽¹⁾ - *Hisashi Nagano*⁽²⁾ - *Hiroki Mizuno*⁽²⁾ - *Takahiro Ito*⁽²⁾ - *Tatsuo Nojiri*⁽²⁾ - *Miyako Iida*⁽¹⁾ - *Hitoshi Maruyama*⁽³⁾

Hitachi, Ltd., Research & Development Group, Tokyo, Japan (1) - Hitachi, Ltd., Service & Platform Business Unit, Ibaraki, Japan (2) - Hitachi, Ltd., Railway Systems Business Unit, Tokyo, Japan (3)

Keywords: Trace explosives detection, mass spectrometry, high-throughput, TNT, TATP

Introduction

The threats of terrorism using explosives have been a serious problem around the world [1]. In this study, we developed two types of mass spectrometry-based explosive trace detectors (ETDs), namely, a high-throughput ETD which detects explosive particles using non-contact sampling, and a compact ETD which detects explosive particles and their vapor by one sampling unit.

Methods

Based on the prototype of high-throughput ETD [2, 3], the non-contact sampling system was improved using a fluid dynamic simulation for stable high sampling efficiency. A mass spectrometer which operated as a sensor in that ETD was downsized, leading to the downsizing of overall system. As the other type of ETD, we combined such mass spectrometer with both a vapor introduction port and swab sampling function.

Results

The newly developed high-throughput ETD has an ID card insertion port where several air nozzles are located. The position of such nozzles were decided based on the fluid dynamic simulation. An air jet from the nozzles efficiently dislodge explosive particles from the ID card and explosive molecules are then detected by mass spectrometry. We confirmed that the high-throughput ETD detected one of representative explosive molecule, TNT within 3 sec after the insertion of the ID card and its sensitivity was ng level. This throughput enables that ETD to operate as a security gate inspecting all people passing the gate. The compact ETD with a weight of around 60 kg has a function of continuous monitoring of surrounding gas and swab sampling. Experimental results suggested that the compact ETD detected vapor of TATP, which is popular improvised explosive device, at a concentration of below 0.1 ppm and TNT particle at ng level. The compact ETD can be interspersed in several locations to increase the security level in a wide area.

Conclusions

High-throughput ETD and compact ETD were developed based on the same mass spectrometer. Use of those ETD as a security gate and an area monitoring system is considered to increase the security level in a wide area of the city.

Novel Aspect

Based on the mass spectrometer as a core technology, changes in a sampling system resulted in two types of ETD used for several purposes.

References

1. United States Department of State Publication, "Country reports on terrorism 2016", <https://www.state.gov/j/ct/rls/crt/2016/index.htm>
2. Hashimoto Y., Nagano H., Takada Y., Kashima H., Sugaya M., Terada K., Sakairi M., Rapid Commun. Mass Spectrom, 28(12), 1376-1380 (2014).
3. Takada Y., Nagano H., Kawaguchi Y., Kashima H., Sugaya M., Terada K., Hashimoto Y., Sakairi M., IEEE Sens. J., 16(5), 1119-1129 (2015).

MP-255 / USING AUTOMATED LOW VOLUME PIPETTING TO ENABLE HIGH THROUGHPUT DART-MS DETECTION OF DRUGS DIRECTLY FROM BIOLOGICAL MATRICES

Frederick Li - Paul Liang - Brittany Laramée - Brian Musselman

IonSense, Inc., Chemistry, SAUGUS, United States

Keywords: DART-MS, Drugs of abuse, Automation, High-throughput Experimentation, Toxicology

Introduction

Direct Analysis in Real Time Mass Spectrometry (DART-MS) is a technique capable of directly analyzing biological samples without sample preparation. Utilizing a 384 format nanoliter pipettor to spot nanoliter droplets on a wire mesh, matrix effects can be reduced up to a factor of 100 times, thus enhancing analyte detection. The 384 sample mesh can then be analyzed with DART-MS at the rate of one sample per second, achieving complete analysis in less than 10 minutes.

Methods

Urine was spiked with various drugs and internal standard (IS) and serially diluted to simulate real-world urine samples submitted for drug tests. Nanoliter droplets of the spiked urine were spotted on a mesh screen using a 384 nanoliter pipettor. A DART ion-source with a high throughput (HT) module was interfaced to a quadrupole mass spectrometer to enable rapid analysis of the mesh screen samples.

Results

Various sample deposition volumes, 50 – 1000 nL, and sample concentrations, 50 ng mL⁻¹ – 1000 ng mL⁻¹, were evaluated to determine the level of matrix effect reduction as well as the lowest droplet volume that provided sufficient response for detection. A constant final amount of material was deposited using different size droplets of varying concentration; Preliminary results showed that the low droplet volume provided comparable signal intensity and a significant reduction of matrix effects. A five-point calibration curve was constructed using 200 nL droplets, and quantification with IS correction showed that the data was reproducible, with a RSD% less than 1%. The drug concentrations were accurately quantified between 50 – 1000 ng mL⁻¹ within a 10% error.

Conclusions

Using nanoliter volume quantities of sample reduced the matrix interference and enabled direct analysis of the samples with DART-MS without any prior sample pre-treatment. The high throughput DART-MS configuration can enable the analysis of 384 samples in approximately 7 minutes and is capable of providing accurate and reproducible quantitation. Future work will include a wider drug panel and different biological matrices.

Novel Aspect

Nanoliter volumes can reduce the matrix effects of biological samples, therefore enabling rapid analysis of biological samples using DART-MS without sample manipulation.

MP-256 / TOWARD HIGH THROUGHPUT SCREENING OF CELL-MATERIAL RESPONSE USING LIQUID EXTRACTION SURFACE ANALYSIS-MASS SPECTROMETRY (LESA-MS)

Joris Meurs⁽¹⁾ - *Josephine Bunch*⁽²⁾ - *Morgan Alexander*⁽¹⁾ - *David Barrett*⁽¹⁾ - *Dong-Hyun Kim*⁽¹⁾

Advanced Materials & Healthcare Technologies Division, University of Nottingham, Nottingham, United Kingdom (1) - National Centre for Excellence in Mass Spectrometry Imaging, National Physical Laboratory, Teddington, United Kingdom (2)

Keywords: Cell-material response, high-throughput screening, in situ trypsin digestion, protein analysis, LESA-MS

Introduction

Pluripotent stem cells are a valuable source for cell production. Biomaterials have been discovered that

assist in pluripotency maintenance, but the response mechanism is not well understood. Pre-adsorption of proteins is crucial in this process [1]. LESA-MS has shown potential to be used for analysis of adsorbed proteins [2]. Here, we describe the development of a high throughput LESA-MS strategy to analyze strongly attached proteins in a complex matrix.

Methods

Essential 8™ medium (E8) was incubated for 1 hour onto the Droplet Microarray which has shown to allow robust automatic surface extraction [3]. Removal of interfering compounds was optimized using ammonium formate solutions and a 10-minute washing. Next, optimization of in situ trypsin digestion was done using organic-aqueous systems and MS-compatible detergents which have shown to enhance in-solution digestion [4]. Tryptic digests were then analyzed by LESA-MS.

Results

Surface washing was optimized by varying the NH_4HCOO concentration. To determine the washing efficiency and reproducibility, the total sodium intensity was calculated using an in-house MATLAB script. One washing round significantly reduced the amount of sodium independent of the NH_4HCOO concentration. More washing rounds as well as the NH_4HCOO concentration were not found to have a significant effect on the washing efficiency. The most reproducible washing was obtained using deionized water. Performing trypsin digestion with a mixed organic-aqueous system (DMF/50 mM NH_4HCO_3 (30:70 v/v) increased the number of identified peptides as well as the MASCOT scores for BSA compared with trypsin digestion using an ammonium bicarbonate solution solely. Addition of a MS-compatible detergent seems to interfere with in situ digestion or MS analysis due to insufficient removal of the breakdown products.

Conclusions

We showed that with discrete sample areas the extraction repeatability and detection limits of LESA-MS for peptides can be improved. Further, using an optimized washing protocol inorganic salts present in E8 could effectively be removed prior to in situ digestion enhancing the peptide signal. In addition, the use of organic-aqueous trypsin solution significantly enhances the in situ digestion efficiency and protein identification.

Novel Aspect

This work is the first to show a strategy directing towards high throughput analysis of surface-adsorbed proteins for media- and cell-derived proteins.

References

1. Celiz A. D., Smith J. G. W., Langer R., Anderson D. G., Winkler D. A., Barrett D. A., Davies M. C., Young L. E., Denning C., & Alexander M. R., *Nature Materials*,13,570-579 (2014).
2. Rao W., Celiz A. D., Scurr D. J., Alexander M. R., & Barrett D. A., *Journal of The American Society for Mass Spectrometry*,24,1927-1936 (2013).
3. Meurs J., Alexander M. R., Levkin P. A., Widmaier S., Bunch J, Barrett D. A., & Kim D. -H., *Analytical Chemistry*,90 (10), 6001–6005 (2018).
4. Chen E., Cociorva D., Norris J. L., Yates J. R. III, *Journal of Proteome Research*,6,2529-2538 (2006).

MP-257 / DEVELOPMENT OF A MULTIMODAL IMAGING PLATFORM BASED ON PHOTOINDUCED THERMAL DESORPTION AND ANALYSIS BY ATMOSPHERIC PRESSURE CHEMICAL IONIZATION MASS SPECTROMETRY

Matthias Lorenz^(1,2) - **Mario Viani**⁽³⁾ - **Aleksander Labuda**⁽³⁾ - **Stephen Jesse**⁽²⁾ - **Roger Proksch**⁽³⁾ - **Olga Ovchinnikova**⁽²⁾

University of Tennessee, Knoxville (1) - Oak Ridge National Laboratory, Center for Nanophase Materials Sciences, Oak Ridge, United States (2) - Oxford Instruments, Asylum Research, Santa Barbara, United States (3)

The key to advancing materials is to understand and control their structure and chemistry. Thorough chemical characterization can be challenging since many existing techniques analyze only a few properties of the specimen, thereby requiring multiple measurement platforms to acquire the necessary information. The multimodal combination of atomic force microscopy (AFM) and mass spectrometry (MS) transcends existing analytical capabilities for nanometer scale spatially resolved correlation of the chemical and physical properties of a sample surface. We recently introduced the utilization of a photoinduced cantilever heating technology developed by Oxford Instruments for the localized thermal desorption and analysis and presented a closed cell design for sampling on an Oxford Instruments Cypher ES microscope to interface with a Thermo Orbitrap Velos Pro mass spectrometer using inline atmospheric pressure chemical ionization (APCI). The photoinduced cantilever heating technology works with standard AFM probes that are compatible with advanced AFM modes for the nanomechanical and electromechanical characterization of samples. We previously demonstrated below 500 nm spatial resolution for the spot-sampling by thermal desorption from thin layers and the chemical analysis of small organic molecules in full scan MS mode.

We demonstrate the application of multiple and advanced AFM measurement modes such as AMFM nanomechanical characterization on a single AFM cantilever, combined with photoinduced thermal desorption and analysis by mass spectrometry to link chemical composition with material functionality. We show the chemical analysis by mass spectrometry of gas phase species evolving from polymeric material in contact with a heated AFM probe, spatially resolved with nanometer resolution, and identify small organic molecules and characteristic fragments and pyrolysis products from the polymer. We present results from systematic studies of factors limiting the efficiency of transport and ionization of material evolving from the sample surface, with the objective to enhance the achievable spatial resolution and compound coverage of the multimodal imaging platform. Parameters studied include the timing of the photothermal cantilever heating and analysis by mass spectrometry, transport conditions for the gaseous material from the closed desorption cell to the mass spectrometer, and the operating conditions of the inline ionization stage. Computational fluid dynamics (CFD) simulations were used to study the uptake and transport of material. Experimental studies on prototype cell models prepared by additive manufacturing provide support for the modeling data.

MP-258 / WHAT DETERMINES THE POST-IONISATION EFFICIENCY IN MALDI-2: A COMBINED SOFTWARE/HARDWARE-BASED SET-UP TO SEMI-AUTOMATICALLY CHARACTERIZE THE ROLE OF FIVE RELEVANT INPUT PARAMETERS

Alexander Potthoff - Klaus Dreisewerd - Jens Soltwisch

Institute of Hygiene, Uni Münster, Münster, Germany

Keywords: MALDI-2, Post-ionisation, Laser fluence, wavelength

Introduction

MALDI with laser-induced post-ionization (PI), also named MALDI-2, was recently introduced by us to increase ion yields for numerous classes of lipids, metabolites and carbohydrates [1]. Here we present a semi-automatic protocol to investigate and optimize the efficiency of the PI process in dependence of five relevant input parameters: pulse energies of the two lasers, PI laser wavelength, delay between the laser pulses, and buffer gas pressure in the ion source.

Methods

All experiments were conducted with a modified Synapt G2-S mass spectrometer (Waters) [1]. A N₂ laser (337 nm) served for desorption/ablation of material and a wavelength-tunable optical parametric oscillator (OPO) laser (relevant range, 220 – 350 nm) for PI. 2,5-Dihydroxybenzoic acid (DHB) and a mixture of polar phospholipids were investigated as matrix/analyte system. Samples were prepared by spraying the prepared mixture onto a glass slide.

Results

At the core of our study stands a newly developed setup that combines hardware components and LabVIEW VI-based software for control and read-out of all relevant input parameters. In this way, during a MALDI imaging run on a homogenous sample, up to four (out of the overall five studied) input parameters were systematically varied and mass spectra recorded simultaneously. For example, by use of motor-driven laser energy attenuators (neutral gradient filters) and solenoid gas valves, laser fluences and buffer gas pressure in the ion source were precisely adjusted in predefined steps. The complex obtained MS data set including total ion counts and abundances of molecular ion signals of matrix, analytes and fragments in dependence of the respective input parameters was used to identify optimal MALDI-2 conditions for the tested DHB/lipid system. As will also be discussed in our poster, the so collected multidimensional data cube provides valuable insight into the mechanisms underlying the MALDI and MALDI-2 processes [2].

Conclusions

Our new software/hardware setup allowed for (rapid) characterisation of the MALDI-2 process in dependence of all relevant input parameters. Our approach could not only support further optimization of MALDI-2 for different matrix systems for lipids, but also allow for a screening of new analytes (e.g., peptides and proteins) which under current conditions show no significant PI. Furthermore the approach is a valuable tool to study the MALDI and PI processes.

Novel Aspect

Combined hard-/software setup for optimization and characterization of the MALDI-2 process in dependence of five input parameters and preliminary data for a DHB/lipid system.

References

1. Soltwisch et al., Science, 348 (6231), 211-215 (2015)
2. Niehaus et al., Sci. Rep., 8, 7755 (2018)

MPS-S10 - Polymers

Chairs: Marek M. Kowalczyk, Leopoldo Ceraulo

MP-259 / A COMBINATION OF FAIMS SEPARATION WITH SYNCHROTRON UV-PHOTO ACTIVATION FOR THE TANDEM MS ANALYSIS OF SYNTHETIC POLYMERS

Véronique Legros⁽¹⁾ - **Ines Aloui**⁽¹⁾ - **Alexandre Giuliani**⁽²⁾ - **William Buchmann**⁽¹⁾

LAMBE (CNRS-UEVE-CEA-UCP, UMR8587), Université Paris-Saclay (Univ Evry), Evry (1) - Synchrotron SOLEIL, UAR1008, CEPIA, INRA, Saint-Aubin (2)

Keywords: synthetic polymers, electrospray, tandem mass spectrometry, UV photo-activation, field asymmetric waveform ion mobility

Introduction

Mass spectrometry (MS) of synthetic polymers most often uses single stage MS alone but characterization of polymers can remain ambiguous without complimentary tandem MS experiments [1,2]. UV-photo-activation creates new opportunities for the characterization of polymers [3,4]. The goal of this work was to evaluate the performances of wavelength-tunable UV photo-activation of synthetic polymers, combined or not with field asymmetric waveform ion mobility.

Methods

Precursor ions produced by ESI (positive mode) from a variety of commercial polymers (polyethylene glycol (PEG), polymethylmetacrylate (PMMA), and polysiloxane (PDMS), Sigma-Aldrich), have been stored in a linear ion trap (LTQ Thermo), and then submitted to UV synchrotron irradiation during a variable time period and over a range of wavelengths. For the highest mass polymers, UV activation was combined with field asymmetric waveform ion mobility spectrometry.

Results

By applying suitable compensation voltages in FAIMS, it was possible to specifically isolate certain $(M+n\text{Cation})^{n+}$ adducts produced by ESI and to avoid the overlapping problems. The selected ions could be trapped in a linear ion trap, and then submitted to UV synchrotron irradiation during a variable time period, at different photon energies. For example, the activation of $[\text{PEG95}+4\text{Na}]^{4+}$, produced different effects as a function of the activation method (low energy CID vs UV-Photoactivation 16eV, 10s). As low energy CID only led to the loss of one sodium cation according to the equation $[\text{PEG95}+4\text{Na}]^{4+} \rightarrow [\text{PEG95}+3\text{Na}]^{3+}$, the UV irradiation $[\text{PEG95}+4\text{Na}]^{4+}$ led to a wide series of fragment ions. Polydimethylsiloxanes also gave net differences. Series of UV-photoactivation mass spectra, recorded with variable photon energies (6 to 24eV), showed that the fragmentation started with photon energies slightly above the supposed ionization potential of the oligomer. Important fragmentations were observed at 16eV.

Conclusions

In this work, we demonstrated the potential of UV photo-activation of ions, produced from polymers and stored in a linear ion trap. The most useful UV wavelengths were identified. The use of FAIMS made it possible the selection of certain precursor ions that were difficult to isolate (particularly useful for highly charged ions). UV photo-activation could lead to different and useful fragmentations compared to the low energy CID of the same precursors.

Novel Aspect

The use of synchrotron UV beam for tandem MS experiments from polymers and the combined use of ion mobility for multiply charged precursor ions selection.

References

1. Wesdemiotis, C. ; Solak, N. ; Polce, M. J. ; Dabney, D. E. ; Chaicharoen, K. ; Katzenmeyer, B. C. *Mass Spectrometry Reviews*, 30, 523(2011).
2. Terrier, P. ; Buchmann, W. ; Desmazieres, B. ; Tortajada, J. *Analytical Chemistry*, 78, 1801(2006).
3. Girod, M. ; Antoine, R. ; Lemoine, J. ; Dugourd, P. ; Charles, L. *Rapid Communications in Mass Spectrometry*, 25, 3259 (2011).
4. Girod, M. ; Antoine, R. ; Lemoine, J. ; Dugourd, P. ; Charles, L. *International Journal of Mass Spectrometry*, 333, 27 (2013).

MP-260 / "ADVANCED KENDRICK MASS DEFECT ANALYSIS": A TOOLKIT

Thierry Fouquet⁽¹⁾ - ***Robert B. Cody***⁽²⁾ - ***Takaya Satoh***⁽³⁾ - ***Hiroaki Sato***⁽¹⁾

National Institute of Advanced Industrial Science and Technology (AIST), Research Institute for Sustainable Chemistry (ISC), Tsukuba, Japan (1) - JEOL USA, Peabody, USA (2) - JEOL, Akishima, Japan (3)

Keywords: Data processing, Kendrick mass defect, resolution-enhanced analysis, high-resolution mass spectrometry, polymers

Introduction

Developed in the eighties from a founding change of mass scale proposed by Kendrick fifty years ago, the

Kendrick mass defect (KMD) analysis has been advantageously used for an easy graphical interpretation of high-resolution mass spectra from carbonaceous samples. Recently extended to the processing of data from polymers [1], the traditional KMD analysis is now enjoying a new lease of life with the introduction of several concepts for an unprecedented visualization experience.

Methods

Regular Kendrick masses (KMs) are $KM(R) = m/z * \text{round}(R)/R$ with m/z the accurate mass-to-charge ratio and R the exact mass of a repeating unit (IUPAC scale). “Advanced” KMs compiling a “resolution-enhancement” [2-4] and a “charge-dependence” [5] are $KM(R,Z,n) = Z * m/z * (\text{round}(R)+n)/R$ with n and Z being positive/negative integers (calculations using spreadsheet or dedicated programs). KMDs and the “remainders of KMs” [6] are $KMD = \text{round}(KM) - KM$ and $RKM(R) = KM / \text{round}(R) - \text{floor}(KM / \text{round}(R))$.

Results

Beyond the simplest KMD plot (KMDs. vs. m/z) from the mass spectrum of a singly charged homopolymer, a regular KMD plot from multiply charged (co)polymer ions displays the charge state distribution instantly via an isotopic split. To overcome possible isotopic misalignments, a new “charge-dependent” KMD plot clusters the ion series at a chosen charge state into a single cloud facilitating their filtering and assignment. The new resolution-enhanced KMD plot strikingly improves the separation of ion series taking full advantage of the spectral width. Using well-fitted parameters, mass spectra of copolymers become a set of elliptic clusters discriminating ions by their end-groups and co-monomeric content. Terpolymers are also easily turned into a series of copolymeric clusters in a first resolution-enhanced plot, further filtered and re-analyzed via a second resolution-enhanced KMD plot or a “degree of polymerization” plot. The new “remainders of KMs” allow mass spectral data of low mass accuracy such as MALDI-TOF/TOF spectra to be turned into informative plots (RKMs. vs. m/z).

Conclusions

The “Advanced KMD analysis” broadens the scope of applications of a very simple data processing tool to low or high mass accuracy data, to singly or multiply charged ions and to homo- / co- / terpolymer ions with one appropriate KMD plot in every case. With a rationalized approach to help users choosing the well-fitted parameters, it may be now routinely used for a graphical and highly visual exploration of mass spectral data.

Novel Aspect

From the historical definition of a Kendrick mass ($KM = m/z * 14/14.0157$) to a resolution-enhanced and charge-dependent KM, we propose a deeply rejuvenated KMD analysis.

References

1. Sato H., Nakamura S., Teramoto K., Sato T., J. Am. Soc. Mass Spectrom. 25, 1346 (2014).
2. Fouquet T., Sato H., Mass Spectrom. (Tokyo), 6, A0055 (2017).
3. Fouquet T., Sato H., Rapid Commun. Mass Spectrom. 31, 1067 (2017).
4. Fouquet T., Sato H., Anal. Chem. 89, 2682 (2017).
5. Fouquet T., Cody R. B., Ozeki Y., Kitagawa S., Ohtani H., Sato H. J. Am. Soc. Mass Spectrom. und rev (2018).
6. Fouquet T., Satoh, T., Sato H., Anal. Chem. 90, 2404 (2018).

MP-261 / POTENTIAL OF LIQUID EXTRACTION SURFACE ANALYSIS FOR MASS SPECTROMETRIC EVALUATION OF SAFETY AND FUNCTIONALIZATION OF POLYMER FOOD PACKING MATERIALS.

Ambre Issart - Simon Godin - Hugues Preud'homme - Ahmed Allal - Joanna Szpunar

IPREM- UMR 5254 CNRS- UPPA – Université de Pau et des Pays de l'Adour – Hélioparc – 2av. du president, Angot 64053 Pau Cedex, France

Keywords: Surface analysis, liquid extraction surface analysis, polymers additives, degradation products, food packaging.

Introduction

The migration of constituents of polymer packaging into food has two aspects. Additives, post-polymerization or degradation residues, Non-Intentionally Added species (NIAS) are a risk for the consumers but, at the same time, “intelligent” packaging releasing natural antioxidants improve food shelf life. Hence, it is essential to develop methods to evaluate the safety and functionalization of polymer food packaging.

Methods

LESA-MS (liquid extraction surface analysis mass spectrometry) combined with nanoESI ionisation was optimized for the direct analysis of polymer packing materials. The method does not require any sample preparation, is fast, eco-friendly and highly reproducible. It appears to have many advantages over other direct analysis techniques and is well adapted to the polymer analysis.

Results

The method was successfully applied to follow the release of oligomers from additive-free low-density polyethylene (LDPE) during recycling. The increasing number of cycles resulted in the gradual appearance of clusters of species with m/z separated by 14 Da (the monomer) thus proving polymer degradation. LESA-MS analysis also allowed the detection of a series of species corresponding to residual polyethylene oxide - a plasticizing additive used for PET synthesis or recycling. The efficiency of the natural antioxidants (vitamin E and C) on the polymer stabilisation was proven by the comparison of the LESA-MS spectra of pure and functionalized LDPE. Moreover, the release of these additives (meant to protect food packed using this “intelligent” material) could be quantified using laboratory prepared solid state standards.

Conclusions

LESA-MS was optimized to study the release of low molecular species from polymer food packing films and the monitoring of degradation products during recycling of additive-free LDPE and migration of antioxidants from functionalized LDPE. The use of solid state standards made possible quantitative monitoring of the antioxidant release and thus quality control of “intelligent” packaging.

Novel Aspect

LESA-MS is a fast and reliable ecological tool with a high potential for the evaluation of safety and functionalization of polymer food packing films.

MP-262 / STUDY ON THE SYNTHESIS AND CHARACTERIZATION OF SURFACE ACTIVITIES OF HYDROPHILIC DERIVATIVES OF β -SITOSTEROL

Jong Choo Lim

Dongguk University, Department of Chemical and Biochemical Engineering, Seoul

Keywords: β -Sitosterol, Hydrophilic derivatives, Poly(ethylene glycol), Surface activity

Introduction

Phytosterol means all the alcohol compounds with steroid moiety found in higher plant life and at least 44 sterols have been identified [1]. The most abundant phytosterol is β -sitosterol (24-ethyl-5-cholestene-3-ol) which has the greatest potential for the production of steroidal drugs, cosmetics for healing effects on damaged skin and anti-inflammatory effects, and food ingredients that function as cholesterol-lowering

agents. However, hydrophobic and lipophobic properties of β -sitosterol have prevented its widespread application.

Methods

A hydrophilic phytosterol derivative (HPS) was prepared by coupling a hydrophilic matrix carboxyethyl- β -sitosterol (CES) to poly(ethylene glycol) (PEG), where CES was prepared by reacting β -sitosterol with succinic anhydride. The chemical structure and the distribution of isomers of HPSs were investigated by H-NMR and Maldi-ToF respectively. Surface activities of HPSs were also measured and discussed as a function of molecular weight of the parent PEG.

Results

DS (degree of substitution) was calculated from $^1\text{H-NMR}$ data and the DS values of HPSs (product having one sterol moiety at one end of PEG) were close to 1.0 and 2.0 for Di (product having two sterol moieties at both ends of PEG). From the Maldi-ToF spectra of PEG, HPS and Di with different chain lengths, each peak has intervals of 44, which is equivalent to the mass of monomeric unit of PEG, and all spectra shows a narrow range in molecular weight. Spectra of HPSs displays identical pattern with those of corresponding PEGs, except that the molecular weight (m/z) of HPS shifts to the higher m/z to the extent of approximately 500 m/z compared to that of PEG. Spectra of Di also displays a displacement to the higher m/z to the extent of approximately 1,000 m/z compared to that of PEG. These results strongly indicate that the HPSs prepared are pure mono-type derivatives. The solubility of HPS in water increased as the molecular weight of the parent PEG increased. However, based on the β -sitosterol moiety in HPS, HPS-1000 prepared from PEG 1000 showed the highest solubility.

Conclusions

Three kinds of HPSs were synthesized from PEG by reacting CES with 5 molar excess of corresponding PEG and removing unreacted PEG by extraction. All HPSs were soluble in water and found to have excellent surface activity and the same cholesterol-lowering effect as β -sitosterol. Thus, the potential application of HPS as a surface active material with biological activities of β -sitosterol such as healing effects on damaged skin or anti-inflammatory effects is expected.

Novel Aspect

New types of HPSs having cholesterol-lowering effects was synthesized by coupling β -sitosterol with a hydrophilic matrix to improve water solubility of β -sitosterol.

References

1. Bean, G. A. Adv. Lipid Res., 11, 193 (1973).

MP-263 / MASS-REMAINDER ANALYSIS (MARA): TO UNDERSTAND MASS SPECTRA OF COPOLYMERS

Tibor Nagy - Ákos Kuki - György Deák - Lajos Nagy - Miklós Zsuga - Sándor Kéki

University of Debrecen, Department of Applied Chemistry, Debrecen

Keywords: Mass Remainder Analysis (MARA), copolymer, compositional drift, end group

Introduction

The characterization and the assignment of individual peaks of copolymers are great challenges in the mass spectrometry owing to the huge number of m/z peaks. The Kendrick mass defect (KMD) analysis [1-2] partly solved this issue with limitations such as the high resolution and high accuracy. However, our developed method that is based on the mass remainder after dividing by the exact mass of one of the repeating unit works with less limitations.

Methods

Measurements were carried out with an Autoflex type MALDI-TOF MS instruments (Bruker Daltonics, Bremen, Germany). DHB (2, 5-dihydroxybenzoic acid) and dithranol (1,8,9-anthracenetriol) were used as matrix, while ionizing agents were sodium trifluoroacetate and silver trifluoroacetate (solvents: MeOH, THF). Ethylene oxide- propylene oxide triblock copolymers were investigated. Thioester type agent was used for the RAFT polymerization of polystyrene [3].

Results

A new data handling method, the Mass-Remainder Analysis (MARA), was developed that use the mass remainder values for the identification of the peaks in copolymer mass spectra and correct the intensities with the defects originated overlaps and isotopic peaks. The mass remainder (MR) values are calculated by division of the measured m/z values by the exact mass of one of the monomer unit. The polymer chains differ only in the number of monomer unit (used for division) have the same MR values. Therefore, these are suitable for grouping, classifying and identifying the copolymers with different composition. For the assignation a defined mass list is used therefore, the deisotoping of the spectrum is also done during the assignation step.

Then in the next step the intensities are corrected. The intensities for the assigned peaks and their isotopologues are calculated and subtracted from the measured ones if it is necessary. Furthermore, the intensities of the identified peaks are corrected with the intensity of corresponding isotopologues.

Conclusions

It was demonstrated that our method was capable of characterizing ethylene oxide- propylene oxide copolymers with higher molecular weight, to determine the number average-, weight average molecular weights, the bivariate distribution and the composition drift [4].

Furthermore, polystyrene was identified to investigate the polymerization mechanism of reversible addition-fragmentation chain-transfer (RAFT) polymerization.

Novel Aspect

The Mass Remainder Analysis (MARA) is a new method which capable of characterizing higher molecular weight copolymers with higher efficiency then previously used methods.

References

1. Kendrick E., Analytical Chemistry, 35, 2146 (1963)
2. Sato H., Nakamura S., Teramoto K., Sato T., Journal of The American Society for Mass Spectrometry, 25, 1346 (2014)
3. Moad G., Chiefari J., Chong Y. K., Kristna J., Mayadunne R. T. A., Postma A., Rizzardo E., Thang S. H., Polymer international, 49, 993 (2000)
4. Montaudo M. S., Macromolecules, 34, 2792 (2001)

Acknowledgement

The work was supported by the GINOP-2.3.2-15-2016-00041 project and by the grant K-116465.

MP-264 / COMPOSITIONAL ANALYSIS OF HIGH MOLECULAR WEIGHT BIODEGRADABLE POLYMERS BY HIGH-RESOLUTION MALDI-TOFMS COMBINING ON-PLATE DEGRADATION WITH RESOLUTION-ENHANCED KMD ANALYSIS

Sayaka Nakamura - Thierry Fouquet - Hiroaki Sato

National Institute of Advanced Industrial Science and Technology, Research Institute for Sustainable Chemistry, Tsukuba

Keywords: on-plate degradation, resolution-enhanced KMD analysis, high-resolution MALDI-TOFMS, copolymer, polymer characterization

Introduction

High-resolution mass spectrometry is notoriously limited to the low mass range (<3 kDa) while industrially important polymers are of high molecular weight (typically >10 kDa). We propose a “on-plate degradation” as an in situ sample preparation to reduce the molecular weight of polyesters using an alkaline solution. The low molecular weight degradation products are then amenable to a mass-analysis taking full advantage of the resolving power of the analyzer.

Methods

Poly(3-hydroxybutyrate-co-3-hydroxyvalerate) (P(3HB-co-3HV)) with >100 kDa has been chosen as a reference aliphatic polyester. A sample solution is deposited on the MALDI target plate followed by a drop of the isotope labeled NaOH-CD₃OD. Removing the excess of alkaline solution with water, the matrix solution is finally added to make co-crystal. MALDI mass spectra are recorded using a high-resolution JMS-S3000 spiral TOF mass spectrometer (JEOL, Japan).

Results

Following the on-plate degradation, the mass analysis of P(3HB-co-3HV) reveals a complex fingerprint in the low mass range (<3 kDa) hence validating the degradation procedure. However, the mass difference between HB / HV repeating units and between carboxyl/methyl ester chain-ends are identical at 14 Da (CH₂ moiety). To distinguish the variation of HB/HV composition from the variation of end-groups, the peaks associated with the methyl ester termination are shifted of 3 Da via the on-plate degradation using the isotope labeled NaOH-CD₃OD solution. Instead of a time-consuming peak-picking and manual assignment of data, ion series are conveniently filtered in a KMD plot computed from the complex mass spectrum using our innovative resolution-enhanced KMD analysis [1]. Peak lists are finally turned into a “degree of polymerization plot” [2] displaying the discrete co-monomeric composition at first glance. The average HV/HB content also computed from the plot is eventually in good agreement with the reference value provided by the supplier.

Novel Aspect

The compositional analysis of a high-molecular weight polyester is made possible by combining in situ isotope labeled degradation with an innovative data processing.

References

1. Fouquet T., Sato H., Mass Spectrom. (Tokyo), 6 (1), A0055 (2017)
2. Fouquet T., Coby R. B., Sato H., J. Mass Spectrom., 52, 618-624 (2017)

MP-265 / COLD PLASMA TREATMENT OF POLYMER SUBSTRATES FOR SALD ICP MS

Jan Preisler⁽¹⁾ - *Marek Stiborek*⁽¹⁾ - *Masoud Shekargoftar*⁽²⁾ - *Viktor Kanicky*⁽¹⁾ - *Jakub Kelar*⁽²⁾

Masaryk University, Department of Chemistry, Brno (1) - Masaryk University, Department of Physical Electronics, Brno (2)

Keywords: SALD ICP MS, dried droplet, polymer treatment, cold plasma, DCSBD

Introduction

Substrate-assisted laser desorption inductively coupled plasma mass spectrometry (SALD ICP MS) is a technique for inorganic analysis of liquid samples of sub-microliter volumes in the form of dried droplets on suitable polymer substrates [1]. Here we study influence of substrate surface aging effects on analysis of

cadmium and indium in samples prepared on polymer substrates treated by diffuse coplanar surface barrier discharge (DCSBD) [2].

Methods

A DCSBD device (RPS400-Roplass plasma system 400 W, Roplass, Czech Republic) was used for treatment of polymer substrates at 400W. Contact angle of the deposited droplet with Cd and In was measured using See System E, Advex Instruments, Czech Republic. Samples were desorbed in ablation cell (UP213, New Wave Research, USA) with a 213-nm laser, and content of ^{111}Cd and ^{115}In analyzed by a quadrupole ICP mass spectrometer (7500ce ICP-MS, Agilent, USA).

Results

Four polymer substrates (polyethylene terephthalate glycol, polyethylene terephthalate, polycarbonate and styrene-acrylonitrile) were exposed to cold plasma for 1, 5 and 10 s. Droplets of Cd and In solutions in citrate, biological culture medium and cell lysate were deposited in five replicates immediately, 1, 3, 6 days, 1, 3 or 6 months after the treatment on the polymer substrates. Contact angles of droplets as well as size of the generated stains were characterized. Typically, a longer treatment of surface substrate resulted in lower contact angle. The stain diameter gradually decreased with time after treatment, sometimes showing a maximum at deposition 1 – 6 days after the treatment. Also, the longer the surface treatment was, the slower return to the original state (high value of contact angle) was observed. Citrate concentrations above 5 mM resulted in significant drop of ^{111}Cd intensity. Using In as an internal standard for determination of Cd alleviated the observed matrix and surface related phenomena and reduced RSD by a factor of ~ 2 .

Conclusions

Aging of surface treated by DCSBD has been monitored using physical methods and related to SALD ICP MS response of samples prepared on the surface. The initial experiments show significant changes in contact angle, stain size and SALD ICP MS response. Using an internal standard, RSD values were reduced to $<10\%$. A half-year monitoring is under progress; a procedure for surface treatment and sample preparation will be recommended after completing the entire set of analysis.

Novel Aspect

Substrate surface treatment by DCSBD for dried droplet preparation. Systematic characterization of surface aging and its consequences for SALD ICP MS reproducibility.

References

1. Peš O., Jungová P., Vyhnánek R., Vaculovič T., Kanický V., Preisler J., Anal. Chem. 80, 8725 (2008).
2. Černák M., Černáková L., Hudec I., Kováček D., Zahoranová A., Eur. Phys. J. -Appl. Phys. 47, 1 (2009).

We acknowledge the financial support of the Czech Science Foundation (GA18-16583S) and the Ministry of Education, Youth and Sports of the Czech Republic under the project CEITEC 2020 (LQ1601).

MP-266 / REMAINDERS OF KENDRICK MASS: KEEPING THE BEST OF A KENDRICK MASS DEFECT ANALYSIS FOR LOW-ACCURACY MASS SPECTRA OF POLYMERS

Hiroaki Sato - Thierry Fouquet

National Institute of Advanced Industrial Science and Technology (AIST), Research Institute for Sustainable Chemistry, Tsukuba

Keywords: low mass accuracy, Kendrick mass defect, Remainers, MALDI-TOF/TOF, polymers

Introduction

The traditional Kendrick mass defect (KMD) analysis [1] is restricted to high-resolution mass spectrometry

(low mass range, high-accuracy analyzers) as it evaluates the difference between rounded and exact m/z of polymer ions. The new resolution-enhanced KMD analysis [2,3] makes data in the high mass range or recorded with low-accuracy analyzers interpretable in a KMD plot but point alignment remains of low quality. The “Remainders of KM” (RKMs) overcome this issue.

Methods

Data were recorded using a JMS-S3000 MALDI mass spectrometer (JEOL, Japan) equipped with a high-resolution spiral TOF, a linear TOF with oligomeric resolution and a TOF/TOF device (unitary resolution). A Kendrick mass (KM) is $KM = m/z * (\text{round}(R) + n) / R$ with R the exact mass of the repeating unit and n an integer (n=0 for a regular analysis, n=+/-1,2,3... for a resolution-enhanced analysis). KMDs and RKMs are $KMD = \text{round}(KM) - KM$ and $RKM(R) = KM / \text{round}(R) - \text{floor}(KM / \text{round}(R))$.

Results

The RKM has been found to consist in the difference between a regular KMD and a resolution-enhanced KMD [4]. In other words, RKMs take advantage of the improved separation of ion series from the resolution-enhanced analysis while errors are cancelled from the subtraction to the regular analysis. Point alignments of homologous species in a new RKM plot (RKMs. vs m/z) are perfectly horizontal with an unrivalled separation power. A RKM plot is compatible with a great variety of mass spectral data, starting from single stage mass spectrometry using a high-accuracy analyzer in the high mass range with isotopic resolution. Despite a decrease of the mass accuracy via a loss of resolving power, tandem mass spectra recorded using a TOF/TOF device with unitary resolution down to oligomerically resolved data recorded with a linear TOF are turned into clear RKM plots separating product ion series and distributions via horizontal alignments, respectively [5]. It greatly facilitates the visualization and assignment as well as any subsequent data processing.

Conclusions

The Remainders of Kendrick masses (RKMs) fit onto any sort of mass spectral data from polymeric samples as long as a minimal oligomeric resolution is reached. RKMs also greatly simplify the data processing of copolymeric data [6] and are totally compatible with multiple charging [4]. Linked to the Kendrick mass defects (KMDs) in a unified theory, they are now part of an “Advanced KMD analysis” available in commercial programs computing mass defect plots.

Novel Aspect

Data from low-accuracy analyzers or in the high mass range are turned into informative compositional maps for an easy separation or assignment of polymeric ion series.

References

1. Sato H., Nakamura S., Teramoto K., Sato T., J. Am. Soc. Mass Spectrom. 25, 1346 (2014).
2. Fouquet T., Sato H., Mass Spectrom. (Tokyo), 6, A0055 (2017).
3. Fouquet T., Sato H., Anal. Chem. 89, 2682 (2017).
4. Fouquet T., Cody R. B., Ozeki Y., Kitagawa S., Ohtai H., Sato H. J. Am. Soc. Mass Spectrom. under review.
5. Fouquet T., Satoh T., Sato H., Anal. Chem. 90, 2404 (2018).
6. Nagy T., Kuki A., Zsuga M., Sandor K., Anal. Chem. 90, 3892 (2018).

MP-267 / ANALYSIS OF DEGRADED POLYBUTYLENETEREPHTHALATE (PBT) PRODUCTS BY THERMAL DESORPTION AND PYROLYSIS COMBINED WITH DART-MS (TDP/DART-MS)

Chikako Takei - Kenichi Yoshizawa

BioChromato, Application support, Fujisawa

Keywords: Polybutylene terephthalate (PBT), evaluation of degradation degree, TDP/DART-MS

Introduction

Polybutylene terephthalate(PBT), a thermoplastic engineering elastomer, is widely used in a variety of industrial fields. Due to its industrial applications, it is important to evaluate its degradation degree to facilitate the product development process. The purpose of this work is to evaluate the degradation degree during the initial degradation stage which hasn't been detected by FT-IR, directly and rapidly.

Methods

A new and degraded product made by PBT were used as the samples. Thermal desorption and pyrolysis combined with DART-MS (TDP/DART-MS) was composed of a compact (Bruker) Q-TOF mass spectrometer equipped with a DART ion source (IonSense) and an ionRocket TDP device (BioChromato). The samples were cut to ~1 square millimeter and put into the sample pot. Mass spectra were measured as the samples were gradient heated from room temperature to 600°C at a rate of 100°C/min.

Results

Regarding from the total ion current gram(TIC), as adding the temperature gradient heating by TDP device, thermal desorption and pyrolysis reactions were detected, as shown by the change in intensity. Mass spectra of samples at ca. 220°C (ca. the PBT melting point) were detected the pyrolysis products, monomer and dimer units of PBT. Comparing the extracted ion current gram (EIC) of the PBT monomers (m/z 221.08, $[C_{12}H_{12}O_4+H]^+$) between new and degraded product, there was a difference in the detected temperature and in the amount of PBT monomer. For the degraded sample, the PBT monomer was detected at a temperature much lower than the new sample. Moreover, for the degraded sample, it was detected at much lower than the pyrolysis temperature of the PBT (ca. 400°C). Thus, the PBT monomer that was detected at much lower than the pyrolysis temperature of PBT was presumed a one of the degradation products.

Conclusions

TDP/DART-MS enables detection the initial degradation state which hasn't been detected by FT-IR. Moreover, TDP/DART-MS can be a useful method to evaluate the degradation degree using PBT monomer as a marker. Thus, this method is expected to be applicable to polymer failure analysis, quality control and so on.

Novel Aspect

TDP/DART-MS is useful failure analysis method for thermoplastic engineering elastomers such as PBT.

MP-268 / EVALUATION OF THERMAL HISTORY OF THERMOPLASTIC RESINS BY THERMAL DESORPTION AND PYROLYSIS COMBINED WITH DART-MS (TDP/DART-MS)

Chikako Takei - Kenichi Yoshizawa

BioChromato, Inc., Application support, Fujisawa

Keywords: Polypropylene, Irgafos168, recycle, degradation, TDP/DART-MS

Introduction

Recycled thermoplastic resins are regarded as important materials for realizing sustainable society. Therefore, it is important to control and grasp the thermal history of polymer recycle, because recycled thermoplastic resins are weakened in physical properties by thermal history added. Thus, the purpose of this work was to evaluate the thermal history, directly and rapidly.

Methods

Polypropylene(PP) virgin pellet and 3× recycled pellet were used. Thermal desorption and pyrolysis(TDP) combined with DART-MS was composed of a compact (Bruker) Q-TOF mass spec. equipped with a DART ion

source (IonSense) and an ionRocket TDP device (BioChromato). The samples were cut to ~1 square millimeter and put into the sample pot. Mass spectra were measured as the samples were gradient heated from room temperature to 600°C at a rate of 100°C/min.

Results

Regarding the total ion current gram of the virgin pellet and 3× recycled pellet from room temperature to 600°C, no significant difference was detected, and similar thermal decomposition patterns were observed. However, regarding the mass spectra, the thermal decomposition reaction of polypropylene seemed to begin at 260°C. An antioxidant additive (Irgafos168) and its oxidized compound (tris(2,4-di-tert-butylphenyl) phosphoric acid) were detected at this point, with differences observed between samples. In comparing the intensity ratios of the oxidized and pure between samples, the ratio of the 3× recycled pellet was higher than that of the virgin pellet. This confirmed that the 3× recycled pellet possessed increased oxidized Irgafos168 because of recycle. In addition, in comparing the extracted ion current grams for Irgafos168 between samples, the 3× recycled pellet showed less than the virgin pellet. This confirmed that Irgafos168 was consumed by recycle.

Conclusions

Using TDP/DART-MS, additives contained in PP and PP polymer matrix were detected independently and directly by gradient thermal heating. TDP/DART-MS can be a useful way to evaluate the thermal history of recycled thermoplastic resins quickly and easily, using the antioxidants contained in the thermoplastic resins as markers. Moreover, this method could be applied to quality control for thermoplastic resin products, such as containers and parts.

Novel Aspect

TDP/DART-MS is useful evaluation method of thermal history for thermoplastic resins using the additive marker contained in them.

MP-269 / MALDI-LID-TOF/TOF ANALYSIS OF STATISTICAL AND DIBLOCK POLYACRYLATE COPOLYMERS

James Town - Glen Jones - David Haddleton

University of Warwick, Department of Chemistry, Coventry, United Kingdom

Keywords: MALDI-LID-ToF/ToF, Polymers, MALDI-ToF, Tandem, Fragmentation

Introduction

MALDI-LID-ToF/ToF mass spectrometry is a technique used to aid in the characterization of proteins (1). While MALDI tandem mass spectrometry techniques have been applied extensively to homopolymer compounds (2), with some application to block copolymers (3). Herein we report on MALDI-LID-ToF/ToF as a powerful tool to give information of monomer distributions, and microstructural information, for both block copolymers and statistical copolymers.

Methods

MALDI-ToF and MALDI-LID-ToF/ToF experiments were performed using an UltrafleXtreme MALDI-ToF/ToF analyser (Bruker, Coventry, UK). The MALDI-ToF experiments were performed in reflectron positive ion mode, with an acceleration voltage of 19 kV. MALDI-LID-ToF/ToF experiments were carried out in reflectron positive mode, with an ion source voltage of 8 kV, and a post fragmentation cell acceleration voltage of 19 kV. Polymers were synthesized using photo-RDRP (4).

Results

MALDI-LID-ToF/ToF results for the homopolymers revealed that the lability of the bromine end group led to a fragmentation pattern heavily dependent upon a radical backbiting process. This in turn led to the MALDI-

LID-ToF/ToF spectra only containing a single dominant fragment structure, which had both the alpha end group and an addition CH₂ vinyl, as well as some smaller secondary fragments from the radical backbiting process. This was seen for both the methyl and ethyl acrylate homopolymers.

It is this simplification of the fragmentation pathway to a single fragment that sequencing the copolymers became much less complicated. Simply observing the overall makeup of the spectra granted simple confirmation of the diblock and statistical copolymer structure, as the diblock retained two distinct regions of monomer separation while the statistical copolymer kept a random structure throughout. Evidence of mixing in the diblock also showed the microstructure resolution the technique can obtain.

Conclusions

MALDI-ToF and MALDI-LID-ToF/ToF can provide simple screening of polymer and copolymer compounds, as the techniques use the same sample spot to produce both results, full analysis can be carried out in under an hour. There are clear differences between the diblock copolymer and the statistical copolymer spectra, and, despite the techniques simplicity, determining the co-monomer arrangement of a copolymer is fully within the reach of MALDI-LID-ToF/ToF.

Novel Aspect

Utilizing a technique for copolymer analysis, previously only used for proteins, for fast screening of microstructure determination.

References

1. D. Suckau, A. Resemann, M. Schuerenberg, P. Hufnagel, J. Franzen, A. Holle, *Analytical and Bioanalytical Chemistry*, 2003, 376, 952-965.
2. W. Chrys, S. Nilüfer, P. M. J., D. D. E., C. Kittisak, K. B. C., *Mass Spectrom. Rev.*, 2011, 30, 523-559.
3. A. Maciejczek, V. Mass, K. Rode, H. Pasch, *Polymer*, 2010, 51, 6140-6150.
4. G. R. Jones, R. Whitfield, A. Anastasaki, D. M. Haddleton, *JACS*, 2016, 138, 7346-7352.

TUESDAY, August 28th, 2018

TPS-S01 - Therapeutic drug monitoring and drug discovery

Chairs: Gérard Hopfgartner, Fabio Garofolo

TP-1 / PHOSPHOPROTEOMICS OF NON-SMALL-CELL LUNG CANCER CELLS TREATED WITH ERLOTINIB REVEALS DRUG-RESISTANT SIGNATURES AND POTENTIAL TARGETS

Jun Adachi - Yuichi Abe - Maiko Nagano - Junko Isoyama - Marina Kishida - Takeshi Tomonaga

National Institutes of Biomedical Innovation, Health and Nutrition, Laboratory of Proteomics for Drug Discovery, Osaka

Keywords: peptide fractionation, StageTip, phosphoproteome, non-small cell lung cancer, kinase-substrate enrichment analysis

Introduction

Mass spectrometry-based proteomics is an indispensable tool for molecular, cellular biology and for the emerging field of clinical medicine. Especially, phosphoproteome is one of key signatures to understand the mode of action and mechanism of drug resistance of kinase inhibitors at the molecular level, pathway level

and system level.

Methods

Recently, we have developed a key technology on peptide separation for rapid and sensitive phosphoproteome analysis [1]. We also developed a sensitive method to identify phospho tyrosine proteome using an immune affinity enrichment method [2]. Combining these technologies, we performed temporal characterization of non-small-cell lung cancer cell lines treated with erlotinib.

Results

We quantified over 12000 phosphorylation sites including 600 phosphorylation sites on tyrosine residue across six cell lines. Especially, phosphotyrosine-proteome can detect bypass kinases which were known to be related to the drug resistance, such as FGFR, AXL, MET. Furthermore, network analysis and kinase-substrate enrichment analysis predicted erlotinib-dependent activation of kinases and signal signatures. We extracted kinases and other enzymes which are up-regulated in resistant cells and selected 46 inhibitors for drug screening. 24 of 46 inhibitors inhibited cell growth of at least one resistant cell line.

Conclusions

Combination of phosphoproteome and phosphotyrosine-proteome profiling is useful to detect markers for drug-efficacy and target candidates for overcoming drug resistance.

Novel Aspect

Novel markers for drug-efficacy and target candidates for overcoming drug resistance were identified by phospho proteome analysis.

References

1. J. Adachi et al., *Anal. Chem.*, 88, 7899–7903 (2016)
2. Y. Abe et al., *J. Proteome Res.*, 16, 1077–1086 (2017)

TP-2 / COMPARISON OF METHIONINE OXIDATION PROFILES OF BIOSIMILAR TUR01 MONOCLONAL ANTIBODY TO ITS REFERENCE UNDER OXIDATIVE STRESS

Ahmet Emin Atik⁽¹⁾ - **Z. Zeynep Yildirim Keleş**⁽¹⁾ - **Gülüpek Güven**⁽²⁾ - **İlgü Ece Gülser**⁽¹⁾ - **Yiğit Erdemgil**⁽¹⁾ - **Deniz Bayçın**⁽¹⁾ - **Özge Can**⁽³⁾ - **R. Serdar Alban**⁽¹⁾

Turgut Pharmaceuticals, Biotechnology Group, İstanbul (1) - Acibadem University, Medical Biotechnology, İstanbul (2) - Acibadem University, Department of Medical Engineering, İstanbul (3)

Keywords: mass spectrometry; surface plasmon resonance; biosimilar monoclonal antibody; oxidative stress; oxidation

Introduction

Mass spectrometry (MS) is now routinely used to monitor post-translation modifications (PTMs) levels in biosimilar monoclonal antibody (mAb) characterization. Among them, methionine (Met) oxidation is accepted as one of the most commonly observed PTM during development, manufacturing, formulation and storage processes. To ensure the safety and efficacy of the drugs, Met oxidation should be tightly controlled during mAb development. It may affect biological function, structure and stability of the protein. In the current study, one of the best-selling therapeutic monoclonal antibody and its biosimilar candidate were subjected to the chemical oxidation via different concentrations of hydrogen peroxide. Under forced oxidative stress, Met oxidation profiles as well as the binding kinetics of FcγRI for reference and biosimilar were compared.

Methods

To determine Met oxidation, UPLC-MS experiments were carried out on a Xevo® G2-XS QToF mass spectrometer (Waters, Milford, MA) equipped with electrospray ionization (ESI). Both biosimilar and reference samples were incubated in the dark with different concentrations of H₂O₂ (0.1%, 0.5%, 1.0%, 3.0%) for 24h at room temperature. After applying oxidative stress, both biosimilar and reference antibodies were denatured, reduced, and alkylated before trypsin digestion at 37°C. The resulting peptide mixture were separated by reversed phase C18 column (Waters Acquity UPLC BEH C18, 100 x 2.1 mm, 1.7 µm, 300 Å) with a 90 minute UPLC run. The separated peptides were identified by online tandem mass spectrometry (MS/MS) Q-ToFMS system. Binding kinetics of 0.1% and 3.0% H₂O₂ treated biosimilar and reference antibodies were measured and compared by surface plasmon resonance (SPR) technique on a Biacore T200 (GE Healthcare).

Results

Both biosimilar and reference comprised of one Met on the light chain (LC-Met4) and four Met residues (HC-Met34, HC-Met83, HC-Met256, and HC-Met432) on the heavy chain. LC-MS/MS experiments were used to identify and quantify Met oxidation in the amino acid sequence. Oxidation of Met residue resulted in +16 Da mass shift compared to the unoxidized counterpart. The oxidized peptides were confirmed through careful interpretation of their MS/MS spectra, *i.e.* 16 Da mass increase in b and/or y ions. Incubation of both products with 0.5% H₂O₂ for 24h at room temperature leads to oxidation of Met256 and Met432 residues completely. Those residues are in the Fc region (CH₂-CH₃ interface) and exposed to the protein surface. It has been shown that they are more prone to oxidation compared to the other Met residues. On the other hand, Met34 is located on the complementarity-determining region (CDR) and it may impact antigen binding kinetics. Our data showed that increasing H₂O₂ concentration has nearly the same impact on Met34 oxidation for both reference and its biosimilar candidate. Moreover, Met4 and Met83 residues have shown nearly 20% oxidation under 3.0% H₂O₂ treatment for two products. In addition, the binding kinetics data will also be discussed for oxidized biosimilar and reference mAb.

Conclusions

The high-level similarities in oxidation profiles of Met4, Met34, Met83, Met256 and Met432 residues have shown that the biosimilar and reference mAb have a similar three-dimensional structure. Additionally, SPR analysis has demonstrated that Met34 oxidation has a similar impact on antigen-binding kinetics for both products.

Novel Aspect

Comparison of Met oxidation profiles of biosimilar and reference monoclonal antibodies under oxidative stress through LC-MS/MS and surface plasmon resonance.

TP-3 / QUANTITATIVE EVALUATION OF ADENOSINE 5'-TETRAPHOSPHATE AND OTHER FIVE ANALYTES RELATED TO NICOTINAMIDE PHOSPHORIBOSYLTRANSFERASE BY LC-MS IN MELANOMA CELLS AND MOUSE PLASMA

Michele Bianchi⁽¹⁾ - *Ambra Grolla*⁽¹⁾ - *Cristina Travelli*⁽¹⁾ - *Silvia Garavaglia*⁽¹⁾ - *Armando Genazzani*⁽¹⁾ - *Giuseppe Orsomando*⁽²⁾ - *Erika Del Grosso*⁽¹⁾

University of Piemonte Orientale, Department of Pharmaceutical Sciences, Novara (1) - Polytechnic University of Marche, Department of Agricultural, Food and Environmental Sciences, Ancona (2)

Keywords: Ap4, NAMPT, B16 cells, mouse plasma

Introduction

Nicotinamide phosphoribosyltransferase (NAMPT) is the rate-limiting enzyme in nicotinamide adenine dinucleotide (NAD) synthesis[1]. Adenosine 5'-tetrphosphate (Ap4) is a nucleotide known to be a

vasoactive purinergic mediator in mammals[2]. However, it has never been reported whether NAMPT can catalyze the synthesis of Ap4. So, the aim of the work was to develop a bioanalytical LC-MS method to quantify Ap4 and NAMPT related analytes in biological samples.

Methods

Three cells lines (B16 WT, B16 FLAG-NAMPT and ShNAMPTLow) and twenty-three plasma samples of healthy male C57BL/6 mice were analyzed. The bioanalytical method developed is a LC-MS method that consists in an isocratic elution on a Phenomenex Luna HILIC column. Each analyte (Ap4, adenosine 5'-diphosphate, adenosine 5'-triphosphate, nicotinamide, nicotinamide mononucleotide and NAD) was monitored in MS/MS, SRM or MRM mode.

Results

Intracellular Ap4 levels were increased more than two times in cells over-expressing NAMPT (B16 FLAG-NAMPT: 1.17 ± 0.20 nmol/mg protein v. B16 WT: 0.56 ± 0.10 nmol/mg protein; $p < 0.05$) and were significantly reduced in cells silenced for the enzyme (ShNAMPTLow: 0.31 ± 0.06 nmol/mg protein; $p < 0.05$).

Moreover, metabolic perturbation conditions modulated Ap4 and nicotinamide mononucleotide synthesis in B16 WT cells. In particular, Ap4 and nicotinamide mononucleotide levels were significantly reduced in low glucose concentration and serum free conditions, respectively.

Finally, Ap4 mice plasma levels showed a good correlation with extracellular NAMPT levels, despite their inter- and intra-subject variability. Not surprisingly, Ap4 also correlated, to a greater extent, with ATP levels measured in parallel [3].

Conclusions

A new LC-MS bioanalytical method was developed to quantify six analytes including Ap4 in B16 cells and in two different biological matrices. Murine cells and plasma with engineered or naturally fluctuating NAMPT levels showed matching Ap4 fluctuations.

Novel Aspect

This study reports that NAMPT regulates in vivo Ap4 levels, conferring an unpredicted signaling role to this enzyme.

References

1. Grolla A. A., Torretta S., Gnemmi I., Amoroso A., Orsomando G., Gatti M., Caldarelli A., LimD., Penengo L., Brunelleschi S., Genazzani A. A., Travelli C., Pigment Cell & Melanoma Research 28, 718-729 (2015).
2. Pintor J., Peláez T., Peral A., The Journal of Pharmacology and Experimental Therapeutics 308, 468-473 (2004)
3. Amici A., Grolla A. A., Del Grosso E., Bellini R., Bianchi M., Travelli C., Garavaglia S., Sorci L., Raffaelli N., Ruggieri S., Genazzani A. A., Orsomando G., Cell Chemical Biology, 24, 553-564 (2017)

TP-4 / ANALYTICAL ASPECTS OF SUNITINIB AND ITS GEOMETRIC ISOMERISM TOWARDS THERAPEUTIC DRUG MONITORING IN CLINICAL ROUTINE

Mauro Buzzo^(1,2) - **Bianca Posocco**⁽²⁾ - **Elena Marangon**⁽²⁾ - **Giuseppe Toffoli**⁽²⁾

Clinical and Experimental Pharmacology unit, Aviano National Cancer Institute, Aviano, Italy (1) - Department of Chemical and Pharmaceutical Sciences, University of Trieste, Trieste, Italy (2)

Keywords: TDM, Sunitinib, isomerism, light, HPLC-MS/MS

Introduction

To perform therapeutic drug monitoring (TDM) of Sunitinib (SUN) and its active metabolite *N*-desethyl SUN, a robust and sensitive LC-MS/MS quantification method is required. However, the geometric

photoisomerization of SUN and its metabolite makes working in the dark highly recommended. Unfortunately, these conditions are unsuitable for clinical application and the aim of this work was to study the isomerism phenomenon in human plasma to overcome this issue.

Methods

The geometric isomerism of SUN and *N*-desethyl SUN has been studied in human plasma. A LC-MS/MS method developed and validated for this study, has been applied for the analysis of these samples: the analytes were extracted from the matrix with a protein precipitation and separated on a Synergi Fusion-RP column (Phenomenex, USA) and the system consisted of a Prominence UFLC XR (Shimadzu, Japan) coupled with an API 4000 QTRAP mass spectrometer (SCIEX, USA).

Results

When the SUN samples were exposed to light, the formation of E-isomer did not reach completion with nearly 55% formed after 30 min. An analogue equilibrium achievement was noticed for *N*-desethyl SUN, but its E-isomer reached a maximum of only 22%. Then, the E- to Z-isomer reconversion kinetics was investigated in dark conditions at different temperatures, different analyte concentrations and in presence of acids. The reconversion rate resulted increased by higher temperatures thus indicating a thermal reversion, not influenced by the analyte concentration. Low pH values seemed to stabilize the E-isomer thus slowing down the reconversion. To decrease the E- to Z-isomer reconversion time, samples were incubated in a heated-water bath at different temperatures (40, 50, 60, 70, 80 and 90°C) just before the analysis, reaching the maximum percentage of the Z-isomer for both SUN and *N*-desethyl SUN in only 5 min with temperature equal or higher than 70°C. Moreover, no degradation has been observed for both the analytes after the heating step.

Conclusions

The developed method does not require the protection from the light during the sample preparation but only a heating step with a heated water bath at 70°C for 5 min. The incubation of the samples allows the rapid reconversion of the E-isomer, formed during the sample preparation, to the desired Z-form. Despite the additional step of 5 min little extends the processing time, it strongly facilitates the handling of the samples making this method more feasible for TDM.

Novel Aspect

For the first time, the study of the geometric isomerism of SUN and *N*-desethyl SUN in human plasma has led to the development of an analytical method to support TDM.

TP-5 / A LC-MS-BASED ASSAY FOR THE DETERMINATION OF ACTIVITY AND INHIBITION OF ANOPHELES GAMBIAE 3-HYDROXYKYNURENINE TRANSAMINASE ENZYME

Rossana Canavesi - Ubaldina Galli - Riccardo Miggiano - Franca Rossi - Erika Del Grosso

Università del Piemonte Orientale, Department of Pharmaceutical Sciences, Novara

Keywords: LC-MS/MS, *Anopheles gambiae*, Ag-HKT, Kynurenine pathway, enzymatic assay.

Introduction

In *Anopheles gambiae*, the most efficient vector of malaria parasite *Plasmodium falciparum*, 3-hydroxykynurenine (3-HK) represents a toxic metabolite. In the adult mosquitoes, the excess of this toxic metabolite is removed by a specific 3-HK transaminase (Ag-HKT) which converts the 3-HK into the more stable xanthurenic acid. Interfering with 3-HK metabolism in *A. gambiae* is a potential approach for the development of transmission-blocking drugs and insecticides.

Methods

LC-MS method was optimized for the evaluation of the kinetic parameters of Ag-HKT and the kinetic characterization of inhibitors of the enzyme. An enzymatic assay was set up by using the recombinant Ag-HKT. The chromatographic separation of analytes was performed in a gradient mode on a Synergi column, using an acidic mobile phase. Mass spectrometric detection was achieved with an ion trap equipped with an ESI source, in positive ionization scan, operating in SRM mode.

Results

The LC-MS method was used for the determination of kinetic constants values for Ag-HKT. The K_m and V_{max} calculated for Ag-HKT resulted in agreement with those reported in the literature [1]. Moreover, this selective method was applied for the evaluation of Ag-HKT inhibition by using a compound whose synthesis was previously described in the literature [2] and a new promising inhibitor (3-INI).

The K_i value of the known inhibitor obtained with the LC-MS method was in agreement with that reported in the literature [2]. On the basis of the results obtained in our experiments a competitive mechanism of inhibition was delineated verifying the accuracy and reliability of the LC-MS assay. The method was applied to the characterization of 3-INI, a new synthetic compound. The determination of its mechanism of action through the analysis of the overlaid Lineweaver-Burk reciprocal plot revealed a competitive mechanism of action. The K_i value was determined suggesting that the new synthetic compound is more potent than the reference compound.

Conclusions

The LC-MS method was developed and its applicability was confirmed by determination of the kinetics parameters of Ag-HKT and the determination of inhibitory potency of a known reference compound. Besides, a new synthetic compound was analyzed resulting in a more potent inhibitor with an improved K_i . New compounds will be synthesized starting from the structure of previously studied inhibitors with the aim of discovering more potent and selective Ag-HKT inhibitors.

Novel Aspect

The novelty is the application of LC-MS method for the determination of mechanism of action and K_i values of new potential Ag-HKT inhibitors.

References

1. Rossi F., Lombardo F., Paglino A., Cassani C., Miglio G., Arcà B., Rizzi M., FEBS Journal, 272, 5653 (2005).
2. Rossi F., Garavaglia S., Giovenzana G. B., Arcà B., Li J., Rizzi M., Proc. Natl. Acad. Sci, 103, 5711 (2006).

TP-6 / RAPID, SENSITIVE FENTANYL DETECTION BY THERMAL DESORPTION PTR-MS AND COLLISION INDUCED DISSOCIATION

Luca Cappellin - Felipe Lopez

Tofwerk AG, R&D, Thun

Keywords: drugs, PTR-MS, thermal desorption, collision induced dissociation

Introduction

Drug detection and identification is a current analytical challenge requiring a robust detector which is sensitive, quantitative, and has sufficient separation power to resolve potential interferences. Drugs such as fentanyl are typically used in combination with common street drugs, however because fentanyl is hundred times more potent than heroin it therefore poses unique dangers. An overdose of fentanyl can be caused even by trace amounts.

Methods

PTR-TOF-MS equipped with the a VOCUS reaction cell was tested for high sensitivity drug detection in real-time. The low vapor pressure of many street drugs requires coupling to a thermal desorption system in a swab like fashion. After instrumental calibration, simulated real samples containing few nanograms of fentanyl were tested. Collision Ion Dissociation was induced by varying the extraction voltage to confirm drug identification at structural level.

Results

Fentanyl detection was achieved in seconds without any pre-separation or sample preparation even in complex mixtures. The soft ionization achieved in the Vocus reactor protonated fentanyl ions. Calculated and measured exact mass and isotopic distribution matches perfectly giving support to the molecular composition assignment determined by exact mass analysis. Confirmation of fentanyl detection was achieved by performing real-time Collision Induced Dissociation using VOCUS PTR-TOF. The collisional energy of protonated fentanyl is varied by changing the extraction voltage after the VOCUS reactor. Higher extraction voltages produce higher collisional energies which in turn increase ion fragmentation. The fragmentation pattern provides information about the molecule structure as in MS/MS techniques. At higher collisional energies the parent ion corresponding to protonated fentanyl ($C_{22}H_{28}N_2O \cdot H^+$) is more fragmented into $C_{13}H_{18}N^+$ and $C_8H_9^+$ ions, matching fentanyl fragments from MS/MS databases.

Conclusions

Vocus PTR-TOF coupled to a thermal desorption system allowed real-time detection, identification, and quantification of fentanyl in a swab without any sample preparation or preconcentration step. Collision Induced Dissociation allowed fentanyl confirmation via matching with MS/MS databases. Single digit nanogram detection limits were achieved. Faster desorption rates for example from a swab desorber would result in even lower detection limits.

Novel Aspect

Real-time drug identification and quantification in swabs was achieved by using PTR-TOF equipped with a VOCUS reaction cell and performing Collision Induced Dissociation.

TP-7 / A VALIDATED AND ROBUST ANALYTICAL METHOD FOR PESTICIDES MEASUREMENT IN CANNABIS BY LIQUID CHROMATOGRAPHY TANDEM MASS SPECTROMETRY

Christophe Clarysse

Perkin Elmer, Chrom and Mass Spec Division, Villebon Sur Yvette

Introduction

Quantitative analysis of pesticide residues in cannabis is of great importance to scientific researchers, scientific instruments manufacturers, cannabis industry stakeholders and government regulatory institutions. Besides the complexity of cannabis matrix, more stringent regulatory action limits for pesticides from certain states within the USA, have furthered analytical challenge to accurately quantify these pesticides. Nonetheless, the need for simple, robust and validated methods for pesticide analysis with good data reproducibility and accuracy are critical to all stakeholders in the cannabis industry. This work shows a validated method for analysis of 70 pesticides, including compounds typically analyzed by gas chromatography, in cannabis flower using LC-MS/MS with ESI source. Limits of quantitation were below all current regulatory action limits for pesticides in cannabis in USA. Good reproducibility and precision were obtained for all analytes with %RSD \leq 23 for all pesticides in the cannabis matrix.

Methods

Acetonitrile cannabis extracts were spiked with 70 pesticides at different concentrations and internal standards added. Subsequently, spiked samples were diluted at different levels (1:1 and 1:5) with 0.1%

formic acid in acetonitrile. All diluted and spiked extracts were directly injected into a PerkinElmer LC coupled with the QSight MS/MS system. Calibration standard curves were generated in both diluted standard solvents and the cannabis matrix. Method reproducibility and precision were determined with 7 replicates of spiked sample injections. Limits of quantitation at two different spike levels were established with at least 3 replicate injections. A-20 min LC-MS/MS was used for pesticides analysis in cannabis on a reverse phase C18 column.

Preliminary Results

Due to the complexity of the cannabis matrix, the reproducibility, precision and stability of the instrument are always of utmost importance for any routine analysis of pesticides in a laboratory. In addition, stakeholders are interested in high sample throughput with reliable and validated methods in order to augment their overall efficiency in the cannabis industry. Chromatographic analysis including column re-conditioning for all 70 pesticides was achieved within 20 min. The instrument sensitivity ranged from 0.5 to 50 ppb for all pesticides, which was well below action limits set by California and Oregon. The limit of quantitation for the pesticides in the cannabis matrix was within 0.001-0.1 µg/g for 1:1 diluted spiked samples and 0.001-0.3 µg/g for 1:5 diluted spiked samples. The reproducibility (% RSD; n=7) for each pesticide in both diluted matrices was ≤ 23 %. The calibration correlation coefficient for each analyte in the cannabis matrix was ≥ 0.98 with linear range over 3 orders of magnitude. Inter-day (n=7) variability of area response was ≤ 6 % without the use of internal standards and < 2 % with internal standards.

Conclusions

A high throughput quantitative analytical method for analysis of pesticides in cannabis was fully validated in accordance with the FDA guidelines for chemical analysis in food using liquid chromatography in combination with tandem mass spectrometry. The limit of quantitation for each pesticide was below the required action limits set by California and Oregon states in the USA.

Novel Aspect

Validated liquid chromatography tandem mass spectrometry method for analysis of 70 pesticide residues, including compounds that are conventionally analyzed on gas chromatography platform, in cannabis with limit of quantitation below the action limits set by Oregon and California states.

TP-8 / DATA INDEPENDENT ACQUISITION MODES FOR IDENTIFICATION, QUANTIFICATION AND MONITORING OF LOW-ABUNDANCE HOST CELL PROTEINS DURING MONOCLONAL ANTIBODY BIOPROCESSING

Catalin Doneanu⁽¹⁾ - ***Alex Xenopoulos***⁽²⁾ - ***Romas Skudas***⁽³⁾ - ***Ying Qing Yu***⁽¹⁾ - ***Asish Chakraborty***⁽¹⁾ - ***Weibin Chen***⁽¹⁾

Waters, Biopharmaceutical Sciences, Milford, MA, USA (1) - EMD Millipore, Bioseparations Department, Bedford, MA, USA (2) - Merck, Biopharm Process Solutions, Darmstadt, Germany (3)

Introduction

Quantification of large molecular weight biotherapeutics in biological matrices relies mostly on enzyme-linked immunosorbent assays (ELISAs) mainly due to several advantages: good sensitivity (1-10 ppm), high-throughput, ease-of-use and low cost per sample. In the case of low-abundance protein impurities (HCPs) present in protein therapeutics (e.g. mAbs), these biological assays typically provide the total HCP concentration (usually expressed in ppm or ng HCP/mg mAb), but ELISAs cannot identify and measure individual HCP contaminants. Several MS-based assays have been recently developed to cover this knowledge gap [1-2]. Here we explored the capabilities of a single dimension chromatographic assay, coupled with mass spectrometry using two different data independent acquisition modes for detection of peptide precursors and their fragments.

Methods

The LC/MS assay described here relies on two data-independent acquisition mode implemented on a quadrupole/time-of-flight (QTOF) mass spectrometer. While in MSE mode all peptide ions produced by the electrospray source are transmitted by the quadrupole analyzer, in SONAR mode the quadrupole slides over the mass range of interest during the time required for recording a single MS-spectrum by the TOF-analyzer. Co-eluting precursor ions with different m/z are separated during the rapid quadrupole scan and their corresponding fragmentation spectra are acquired using an identical quadrupole separation. Mass spectra are recorded with high-resolution ($\sim 25,000$) for precursors/fragments.

A monoclonal antibody (mAb), purified by Protein-A affinity-chromatography and then by SCX-chromatography (strong-cation-exchange) using four different elution protocols, was denatured/reduced/alkylated and digested with trypsin.

Preliminary data

LC/MS-assays have recently been adopted as orthogonal techniques to ELISAs for HCP analysis due to their flexibility and sensitivity. Here we describe an efficient analytical scale two step LC/MS assay that allows the identification and quantification of HCPs during mAb purification in a CHO (chinese hamster ovary) cell line. In the first step, called the Discovery Assay, the mAb digest (100-500 μg) is analyzed in SONAR mode on a C18 column (2.1 x 150 mm, operated at 200 $\mu\text{L}/\text{min}$ flow rate) using a 90 min gradient. Two HCPs and 4 spiked protein standards were identified across 5 mAb preparations (one Protein A and 4 SCX fractions) following a CHO database search using Progenesis QIP 4.0 software. The same software was also used to build a spectral library containing 13 CHO HCP peptides and 57 peptides from the spiked protein digestion standards. The spectral library contained the peptide RT, the precursor m/z and the MS/MS fragmentation spectrum of each peptide, thus facilitating subsequent HCP identification and quantification in other mAb preparations.

In the second step of the assay, the same 5 digests were analyzed in triplicate in a different acquisition mode (MSE) using higher throughput runs (30 min gradients, 45 min runs) and lower sample loads (10-20 μg of mAb digest). The MSE acquisition is a data-independent acquisition mode which does not use precursor ion isolation. Instead, in MSE low-energy scans (for precursors) are alternated with high-energy scans (for fragments). The MSE data was searched against the spectral library (containing 70 peptides) to quantify and monitor all 6 previously identified proteins (2 HCPs and 4 spiked proteins). The sensitivity of both assays (discovery and monitoring) was 10 ppm.

SONAR acquisition offers additional selectivity to MSE acquisition, by producing fragmentation spectra with minimized background interferences.

References

1. Doneanu C.E., Xenopoulos A., Fadgen K., Murphy J., Skilton S.J., Prentice H., Stapels M., Chen W., *Mabs* 2012, 4, 24.
2. Doneanu C.E., Anderson M., Williams B.J., Lauber M.A., Chakraborty A., Chen W., *Anal. Chem.* 2015, 87, 10283.

TP-9 / CLINICAL VALIDATION OF A UPLC/MS/MS-BASED ASSAY TO CONCURRENTLY MONITOR FOUR DIRECT ORAL ANTICOAGULANTS IN DRIED BLOOD SPOTS

Kathrin Isabelle Foerster⁽¹⁾ - ***Andrea Huppertz***⁽¹⁾ - ***Andreas D. Meid***⁽¹⁾ - ***Oliver J. Müller***⁽²⁾ - ***Lisa Tilemann***⁽²⁾
- ***Timolaos Rizos***⁽³⁾ - ***Walter E. Haefeli***⁽¹⁾ - ***Jürgen Burhenne***⁽¹⁾

Heidelberg University Hospital, Department of Clinical Pharmacology and Pharmacoepidemiology, Heidelberg (1) - Heidelberg University Hospital, Department of Cardiology, Angiology and Pneumology, Heidelberg (2) - Heidelberg University Hospital, Department of Neurology, Heidelberg (3)

Keywords: Clinical validation, direct oral anticoagulant, dried blood spot, liquid chromatography, tandem

mass spectrometry

Introduction

Monitoring blood concentrations of direct oral anticoagulants (DOAC) is not required regularly [1]. However, concentrations of some DOAC correlate well with (un)desired outcomes [2,3]. Thus, analyzing individual DOAC exposure could be useful for risk assessment, *e.g.* when clearance pathways are impaired. Dried blood spotting (DBS) is a practicable sampling technique for regular monitoring, enabling blood sampling at home at any relevant time point [4].

Methods

DOAC concentrations were quantified using a previously published ultrahigh performance liquid chromatography/tandem mass spectrometry (UPLC/MS/MS) multi-compound assay [5], which was refined for DBS analysis. In a clinical study, approved by the responsible ethics committee, DBS from capillary whole blood and plasma were simultaneously sampled in 33 DOAC patients. Deming regression and Bland Altman analyses were used for evaluation.

Results

Our DBS assay fulfilled pertinent validation standards for a hematocrit range of 33-65 % and was linear within a calibration range of 2.5 (apixaban, rivaroxaban), 4.4 (dabigatran), and 9.3 ng/mL (edoxaban) to 750 ng/mL with only minor ion suppression (matrix effect $\leq 12.5\%$). Inter-assay and intra-assay accuracies ranged between 88.3-110 % and corresponding precisions were $\leq 13.3\%$. DOAC in DBS were stable up to 52 d at room temperature, if stored in the dark and with humidity protection (desiccant). Deming regressions yielded conversion factors for each DOAC (apixaban 1.46, dabigatran 1.51, edoxaban 0.95, and rivaroxaban 1.37), which were used to predict plasma concentrations from (capillary whole blood) DBS. Bland Altman analyses confirmed that plasma concentrations predicted from DBS agreed well with actually analyzed plasma concentrations (71 % for apixaban, 94 % for dabigatran, 100 % for edoxaban, and 88 % for rivaroxaban). Thus, meeting the EMA guideline's demand that in $> 67\%$ of all sample pairs' concentrations must not deviate $> 20\%$.

Conclusions

UPLC/MS/MS for DOAC quantification from DBS constitutes a practicable and sensitive technique to monitor DOAC exposure. The results of the clinical validation study confirmed that blood concentrations from DBS can accurately predict plasma exposure. Thus, the present DBS assay is an attractive option to monitor ambulatory DOAC patients.

Novel Aspect

We present a DBS assay capable of analyzing all DOAC at once, with minor setting requirements regarding sample collection.

References

1. Kirchhof P., Benussi S., Kotecha D., Ahlsson A., Atar D., Casadei B., Castella M., Diener H. C., Heidbuchel H., Hendriks J., Hindricks G., Manolis A. S., Oldgren J., Popescu B. A., Schotten U., Van Putte B., Vardas P., Eur Heart J 37, 2893-2962 (2016)
2. Reilly P. A., Lehr T., Haertter S., Connolly S. J., Yusuf S., Eikelboom J. W., Ezekowitz M. D., Nehmiz G., Wang S., Wallentin L., J Am Coll Cardiol 63, 321-328 (2014)
3. Ruff C. T., Giugliano R. P. ; Braunwald E., Morrow D. A., Murphy S. A., Kuder J. F., Deenadayalu N., Jarolim P., Betcher J., Shi M., Brown K., Patel I., Mercuri M., Antman E. M., Lancet 385, 2288-2295 (2015)
4. Enderle Y., Foerster K., Burhenne J., J Pharm Biomed Anal 130, 231-243 (2016)
5. Foerster K. I., Huppertz A., Müller O. J., Rizos T., Tilemann L., Haefeli W. E., Burhenne J., J Pharm Biomed Anal 148, 238-244 (2018)

TP-10 / DEVELOPMENT OF A SPME TOOL FOR THE DETERMINATION OF A SMALL MOLECULE-DRUG CONJUGATE AGAINST CARBONIC ANHYDRASE IN CANCER CHEMOTHERAPY

Sahar Ghiasikhou⁽¹⁾ - *Jörg Scheuermann*⁽²⁾ - *Samuele Cazzamalli*⁽²⁾ - *Dario Neri*⁽²⁾ - *Renato Zenobi*⁽¹⁾

Department of Chemistry and Applied Bioscience, ETH Zürich, CH-8093 Zürich, Switzerland (1) - Institute of Pharmaceutical Sciences, ETH Zürich, CH-8093 Zürich, Switzerland (2)

Keywords: Capillary gap sampler, ESI Mass spectrometry, Solid phase microextraction, Small molecule drug conjugate, Cancer chemotherapy.

Introduction

The targeted delivery of cytotoxic agents into tissues, especially for malignant cells is an attractive strategy to avoid dose limiting toxicity. Targeting moieties can be antibodies, aptamers and low molecular weight non-peptidic ligands. Previous studies showed the advantages of conjugating small molecules to the drug over attaching macromolecules. One of the common used target moiety for selective delivery of cytotoxic agent is acetazolamide. In this area, there is a pressing need to develop rapid techniques to quantify the amount of drug in the tissue in order to investigate the tumor targeting performance of the different ligands. So-called the “capillary gap sampler” is capable of automated and site-specific extraction of very small sample amounts, which is interesting for this application.

Methods

The capillary gap sampler is a miniaturized sampling device enabling automated low-volume sample handling, for direct interfacing to ESI mass spectrometry. It consists of a liquid bridge of several nanoliters formed between two capillaries, where one acts as the ESI-MS spray needle. Sample extraction is performed by a coated stainless-steel pin which is controlled by a robot arm. Selective drug extraction is possible through immobilizing carbonic anhydrase (CA) on the pin. This modification is performed by overnight reaction between protein with epoxy modified beads which are glued to the tip of the pin. Using fluorescence microscopy, bright field and blue channel images of beads carrying protein linked to 8-anilino-1-naphthalene sulfonate and controls confirmed the attachment of the protein.

Results and discussion

The development started with optimization of the coating procedure, desorption solution, etc. In order to find the optimum desorption solution two tests were performed: in the first one, influence of pH on the binding of acetazolamide with CA was investigated. In the second test, different phases (ACN, MeOH, EtOH, Acetone 60% in water) were used for desorption. ACN was found to be the optimum desorption solution. Extraction of the drug from PBS was performed by dipping the CA modified extraction tool inside the solution. After a quick washing step with water, it enters into the liquid bridge, where the analyte desorbs and is sprayed in the ESI source. Peaks corresponding to the drug (drug^+), ($\text{drug}+\text{Na}^++\text{H}^+$), ($\text{drug}+2\text{H}^+$), (therapeutic warhead=cytotoxic agent coupled to the targeting agent $+\text{H}^+$) were observed in the spectrum. Finally, repeatability of the acetazolamide extraction using CA modified beads are evaluated by performing thirteen 5-minute extractions from 500nM acetazolamide in PBS solution. The relative standard deviation was below 10%, which confirms the repeatability of the method. Drug extraction studies from human plasma will also be shown.

Novel Aspect

Development of an extraction method for evaluation of targeted drug delivery.

TP-11 / IN VIVO ASSESSMENT OF THE CONTRIBUTION OF SURFACTANT REPLACEMENT THERAPY IN RABBITS BY STABLE ISOTOPES NATURAL ABUNDANCE APPROACH

Sonia Giambelluca⁽¹⁾ - *Francesca Ricci*⁽²⁾ - *Manuela Simonato*⁽³⁾ - *Luca Vedovelli*⁽⁴⁾ - *Alessio Correani*⁽⁵⁾ - *Costanza Casiraghi*⁽²⁾ - *Matteo Storti*⁽²⁾ - *Paola Cogo*⁽³⁾ - *Fabrizio Salomone*⁽²⁾ - *Virgilio Paolo Carnielli*⁽⁵⁾

Padova University Hospital, Department of Women's and Children's Health, Padova (1) - Chiesi Farmaceutici, R&D Department, Parma (2) - University of Udine, Department of Medicine, Udine (3) - Fondazione Istituto di Ricerca Pediatrica Città della Speranza, PCare Laboratory, Padova (4) - Polytechnic University of Marche and Azienda Ospedaliero-Universitaria Ospedali Riuniti, Division of Neonatology, Department of Clinical Sciences, Ancona (5)

Keywords: GC-C-IRMS; stable isotope; ^{13}C natural abundance; surfactant; Respiratory Distress Syndrome

Introduction

Respiratory Distress Syndrome (RDS) is a breathing disorder affecting preterms due to lack of lung surfactant, a lipoprotein layer which prevents alveolar collapse [1]. Treatment with exogenous surfactant improves gas exchange and survival [2]. The aim of this study was to estimate the contribution of exogenous surfactant phospholipids to the rabbit alveolar pool in vivo by ^{13}C natural abundance.

Methods

In a feasibility study we measured the $^{13}\text{C}/^{12}\text{C}$ ratio of disaturated-phosphatidylcholine palmitate (DSPC-PA) in poractant alfa and in bronchoalveolar lavages of 20 rabbits. The contribution of exogenous surfactant was then estimated in a rabbit model of RDS [3] including 7 controls and 15 rabbits treated with 50, 100 or 200 mg/kg of poractant alfa. DSPC-PA $^{13}\text{C}/^{12}\text{C}$ ratio was measured by GC-C-IRMS.

Results

The DSPC-PA $^{13}\text{C}/^{12}\text{C}$ ratio of poractant alfa was constant among 7 batches produced over a year with a mean value of -18.8 ± 0.1 ‰ [Range: -18.9 ‰; -18.6 ‰]. The mean $^{13}\text{C}/^{12}\text{C}$ ratio of surfactant DSPC-PA recovered from the lung lavages of 20 rabbits was -28.8 ± 1.2 ‰ [Range: -31.7 ‰; -25.7 ‰]. DSPC-PA in rabbits had a significant lower ^{13}C abundance than poractant alfa, with a mean difference of 10.0‰, wide enough to distinguish the endogenous from the exogenous compound. The contribution of exogenous surfactant to the total alveolar surfactant was calculated in all treated rabbits and it ranged from a minimum of 42.9% in 50 mg/kg group to a maximum of 90.2% in 200 mg/kg group, with a significant dose-effect.

Conclusions

The DSPC-PA of exogenous and endogenous surfactant have a different natural abundance of ^{13}C and this difference could be exploited to estimate in vivo the alveolar contribution of exogenous surfactant in treated rabbits. Since the use of chemically synthesized tracers is not required, it could be useful in human research and in surfactant replacement studies in preterm newborns.

Novel Aspect

This is a novel method to measure in vivo the contribution of exogenous surfactant to the alveolar pool by using stable isotopes at natural abundance.

References

1. Walther F. J., Hernández-Juviel J. M., Waring A. J., Peer J, 2, e403 (2014)
2. Mercier C. E., Soll R. F., Clinics in Perinatology, 20, 711 (1993)
3. Ricci F., Catozzi C., Murgia X., Rosa B., Amidani D., Lorenzini L., Bianco F., Rivetti C., Catinella S., Villetti G., Civelli M., Pioselli B., Dani C., Salomone F., PLoSOne, 12, 1(2017)

TP-12 / BIOPHARMAVIEW™ FOR FAST AND EFFICIENT MONITORING OF pH INDUCED DEAMIDATION

Ferran Sanchez⁽¹⁾ - **Luo Ji**⁽²⁾ - **Kerstin Pohl**⁽³⁾ - **Jörg Dojahn**⁽³⁾ - **Annu Uppal**⁽⁴⁾ - **Sibylle Heidelberger**⁽⁵⁾

SCIEX, Biologics, Barcelona (1) - SCIEX, Biologics, Shanghai (2) - SCIEX, Biologics, Darmstadt (3) - SCIEX, Biologics, Gurugram (4) - SCIEX, Biologics, Warrington (5)

Keywords: deamidation, trastuzumab, attribute monitoring, pH

Introduction

Deamidation is one of the most common modifications which can be induced in monoclonal antibodies during production and prolonged storage. It is critical to evaluate and monitor these biologically relevant product quality attributes (PQAs) as deamidation can cause structural changes, reducing biological activity and efficacy of the biotherapeutic. Here, we monitored the degree and site of deamidation induced in pH stressed antibody digest.

Methods

The trastuzumab samples were exposed to pH stress by increasing the pH of digestion buffer from nearly acidic (6.8) to alkaline (8.8). All measurements were carried out in replicates on using high flow liquid chromatography and a high resolution mass spectrometer using a data dependent (DDA) acquisition strategy. The complete data processing was performed with BioPharmaView™ software using the MAM workflow and reference ranges were set for pass/fail criteria.

Results

The monoclonal antibody digest was exposed to higher pH levels in order to induce deamidation. All asparagine sites were monitored and defined for quantification using MAM workflow in BioPharmaView™ software. The lowest levels of deamidation were observed with the neutral or nearly acidic pH (6.8) as compared to the alkaline pH (8.8) values. The MAM workflow in the BioPharmaView™ software enabled the accurate and automated quantitation of all deamidation sites within the defined criteria of pass/fail using neutral pH values as reference points. The accuracy of the method was assessed by spiking the forced deamidated samples into the control sample at different levels. The sample with the lowest deamidation levels (pH 6.8) was used as control due to lack of samples with nodeamidation. Excellent linearity was observed with the increasing levels of spiking thus verifying the accuracy of the method.

Conclusions

Using BioPharmaView™ software with the attribute calculator for the MAM workflow, levels of deamidation were calculated for all possible deamidation sites within trastuzumab. The data was automatically processed and showed the susceptibility of some residues to pH whereas others are far less sensitive. The results can be used to refine monitoring of attributes to those that are more susceptible to environmental changes.

Novel Aspect

MAM workflow in BioPharmaView™ software for the quantitation of pH induced deamidation.

TP-13 / UNREVEALING PHASE I METABOLISM OF COMBRETASTATIN A4 BY A MULTI-TOOL MS-BASED APPROACH

Karol Jaroach⁽¹⁾ - ***Paulina Zofia Goryńska***⁽¹⁾ - ***Krzysztof Goryński***⁽¹⁾ - ***Tomasz Stefański***⁽²⁾ - ***Barbara Bojko***⁽¹⁾

Nicolaus Copernicus University in Torun, Collegium Medicum in Bydgoszcz, Department of Pharmacodynamics and Molecular Pharmacology, Faculty of Pharmacy, Bydgoszcz (1) - Adam Mickiewicz University, Department of Crystallography, Faculty of Chemistry, Poznań (2)

Keywords: RLM, SPME, combretastatin A4, electrochemistry

Introduction

Drug metabolism remains the main point of interest for clinicians, as comprehensive information regarding metabolic pathways and the mechanisms of action of new drugs are essential for recognition of their possible toxicity and biological activity in patients. The multi-tool analytical approach based on electrochemistry, in silico predictions and in vitro studies with the use of rat liver microsomes were employed for recognition of combretastatin A4 metabolism.

Methods

In silico prediction was carried out using MetaPrint2D-React for recognition of transformations at primary sites of metabolism. For electrochemical reaction the ROXY™ system was used. It conducts redox reactions which are the most common during I phase metabolism. The rat liver microsomes was analyzed via two different extraction approaches protein precipitation and solid phase microextraction. OrbitrapHRMS was used for instrumental analysis.

Results

The use of a high resolution mass spectrometer enabled a non-targeted investigation of metabolic transformations of the parent drug and identification of previously reported [1] as well as unreported metabolites. In total fourteen putative metabolites of combretastatin A4 were found. The main paths of CA4 metabolism are O-demethylation and aromatic hydroxylation and those metabolites were found by all aforementioned metabolism generating methods. There was no differences found between PP and SPME analysis of RLM in terms of identified metabolites. The use of SPME allowed performing time course analysis without need of sample multiplication. The EC reaction resulted in generation 3 metabolites reported also by RLM and one exclusive metabolite not found by other techniques.

Conclusions

The high mass accuracy and fragmentation patterns supported by in silico data showcase the applicability of this approach as an ultimate tool for further identification of putative metabolites detected in microsomes extracts or via EC reactions, especially in cases where authentic standards are not available to confirm the tentative results attained via untargeted analyses.

Novel Aspect

In total of eight metabolites of combretastatin A4 not reported elsewhere were identified via non in vivo based approaches.

References

Aprile, S., Del Grosso, E., Tron, G. C., Grosa, G. (2007). Drug Metabolism and Disposition, 35, 2252-2261 (2007)

Work has been supported by National Science Center No2016/23/N/NZ7//01929

TP-14 / APPLICATION OF TRIPLE QUADRUPOLE MS FOR DRUG-TO-ANTIBODY RATIO AND MOLECULAR WEIGHT DETERMINATION OF ANTIBODY-DRUG CONJUGATES.

Malin Källsten⁽¹⁾ - ***Rafael Hartmann***⁽²⁾ - ***Matthijs Pijnappel***⁽³⁾ - ***Sara Bergström Lind***⁽¹⁾ - ***Jonas Bergquist***⁽¹⁾
- ***Fredrik Lehmann***⁽⁴⁾

Uppsala University, Department of Chemistry-BMC, Uppsala (1) - Uppsala University, Department of Medicinal Chemistry, Uppsala (2) - Recipharm OT Chemistry AB, Uppsala, Sweden (3) - Oncopeptides AB, Stockholm, Sweden (4)

Keywords: Antibody drug conjugate (ADC), RPLC, drug-to-antibody ratio (DAR), triple quadrupole MS, Middle-down

Introduction

Antibody-drug conjugates (ADCs) are an increasingly popular class of therapeutic drugs with an inherently high heterogeneity from a large arsenal of drugs, linkers and antibodies. RPLC-MS using QTOF or Orbitrap analyzers have been used successfully to determine the average drug-to-antibody ratio (DAR) for several constructs. Here we have investigated the potential for the very common triple quadrupole (QqQ) analyzers as an alternative for this application.

Methods

Deglycosylated Trastuzumab and in-house synthesized ADCs generated from Trastuzumab by disulfide reduction with TCEP, followed by conjugation of vcMMAE to the resulting thiols were used for the comparison. Waters Premiere QTOF and Waters XevoTQ-S Micro in combination with Waters Acquity UPLC system and Waters Xbridge C4 column was used for the RPLC separation for the comparison. Mobile phases consisted of 0.5% FA with H₂O or ACN.

Results

Both ADCs and monoclonal antibody (mAb) were analyzed with RPLC-MS on a QTOF and QqQ instruments. Having reduced the disulfide bonds to separate the light- and heavy chains the charge envelopes moves to lower m/z so the loss of ions by a limited m/z range was small. Consequentially the loss of any ions above 2000 m/z in the QqQ was found to have no evident effect on the assigned molecular weights after deconvolution. Comparing the two instruments, preliminary data indicates that the differences in spectral quality related to differences due to age and type of ionization source had bigger effect than the loss of part of the charge envelope on the DAR determination and LOD.

Conclusions

Preliminary data implicates that for average DAR determination of ADCs the QqQ instruments can be a cheaper and more readily available option for smaller labs in early stages of ADC development. The limited m/z range of the QqQ has this far proven to have no major effect on DAR value determination.

Novel Aspect

We have shown that triple quadrupoles can be used for accurate DAR determination of disulfide reduced ADCs and mAbs.

TP-15 / DRAWING THE ACTIVITY OF AURANOFIN BY AN ELECTROCHEMICAL REACTION CELL COUPLED TO AN ESI MS/MS

Pawlak Katarzyna⁽¹⁾ - ***Monika Kupiec***⁽¹⁾ - ***Robert Ziolkowski***⁽²⁾

Warsaw University of Technology, Faculty of Chemistry, Chair of Analytical Chemistry, Warsaw (1) - Warsaw University of Technology, Faculty of Chemistry, Chair of Medical Biotechnology, Warsaw (2)

Keywords: electrochemical chamber, auranofin, redox reaction, pharmacokinetics

Introduction

Gold complexes belong to a group of rapidly emerging novel therapeutic and diagnostic compounds with improved specificity for tumor tissue and strong antiproliferative potency. They are usually responsible for the cell growth-inhibiting effects related to the mitochondrial activation of apoptosis. Auranofin used in the treatment of rheumatoid arthritis was selected for the experiments because of documented activity against ovarian cancer, HIV, E. Histolytica and M. Tuberculosis.

Methods

Redox activity of selected gold(I) complex was investigated by means of a ROXY ERC system (Antec, The Netherlands) consisting of a potentiostat, equipped with an electrochemical micro-reactor cell (μ -PrepCell, Antec, The Netherlands) and an infusion pump (KD Scientific, model KDS100, USA). The Roxy system was online hyphenated to an Agilent 6460 LC/MS ESI triple quadrupole mass spectrometer (Santa Clara, USA) as a detector.

Results

Reduction and intermediate electrochemical conditions corresponding to conditions of cytosol and blood were found optimal to obtain most probable reactive forms for auranofin: (tEP)₂-Au and tetra(tri and di)AtgS-Au, which can interact with bioligands in human organism. Obtained products are in agreement with forms of auranofin bound to GSH and human albumin. Products of deacetylation were minor and observed mainly in the presence of sodium cations. Moreover, sodium ions were found to enhance creation of disulfide bridges during anodic process. ERC-ESI MS results additionally allow to state that deacetylation and loss of Atg or tEP groups by auranofin are irreversible as is disulfide formation. Nitrile group attaching to gold and its influence on transformation of auranofin can mimic the behavior of ligands typical for physiological conditions. Ligand exchange can be proposed as the most probable mechanism for metabolic changes of auranofin.

Conclusions

Electrochemical reaction chamber (ERC) coupled to ESI MS was found to provide significant information considering the transformations of auranofin. Moreover, molecule-specific mass spectrometry delivers complementary data to square-wave voltammetry, which aids in understanding the nature of electrochemical conversions of complex or unstable compounds.

Novel Aspect

For the first time ERC-ESI MS was used to study the transformations of cytotoxic metal transition complexes such as auranofin.

This work was financially supported by the Polish National Center of Science (NCN) within the project DEC-2013/09/B/ST4/00961

TP-16 / THERAPEUTIC EFFECT OF A NATURAL ALKOLOID OBSERVED IN 5XFAD TRANSGENIC ALZHEIMER'S TRANSMODEL: A PROTEOMICS STUDY

İrem Kırış⁽¹⁾ - *Merve Karayel Başar*⁽¹⁾ - *Büşra Gürel*⁽¹⁾ - *Tomasz Mroczek*⁽²⁾ - *Krystyna Skalicka-Woźniak*⁽²⁾ - *Ahmet Tarık Baykal*⁽¹⁾

Acibadem Mehmet Ali Aydınlar University, Medical Biochemistry, Istanbul, Turkey (1) - Lublin University, Pharmacognosy, Lublin, Poland (2)

Keywords: Alzheimer's Disease, neurodegeneration, LC-MS/MS, natural compounds, proteomics

Introduction

Alzheimer's disease (AD), the most common form of dementia (%70). Current drug treatments aim to prevent cognitive decline. No new drug have been approved by FDA since 2003. Studies have revealed that alkaloids with anticholinesterase activity are a promising resource in drug discovery. [1] The aim of the project is to discover natural compounds for the treatment of AD and investigate their effects on neurodegeneration in a comprehensive manner.

Methods

Natural compounds (2 different doses) were administered for 7 days, to 12 m old (pathological phase)

5XFAD transgenic mice (n=8). Littermate mice were used as control group. FDA approved AD drug Galantamine is also used as positive control. Morris water maze (MWM) was performed to investigate changes in learning and memory impairment. Action mechanism of drugs were investigated in 3 different regions of the brain which are cerebellum, cortex and hippocampus via LC-MS/MS based protein expression analysis.

Results and Conclusions

2 of 6 compounds showed significant improvement in learning and memory according to MWM test. The compound which improves the condition most, was used for further proteomic investigation. Protein expression changes between control and drug-administered groups were determined by LC-MS/MS analysis. In each samples, approximately 1250 proteins were identified and differentially expressed proteins between groups were subjected to pathway analysis to enlighten the molecular mechanisms affected by the drug.

Novel Aspect

This is the first study that nLC-MS/MS based proteomics approaches were used in the research of molecular effect mechanism of novel compounds on 5xFAD mouse model.

References

1. Konrath, E. L., Passos, C. d., Klein-Júnior, L. C., Henriques, A. T. J. *Pharm. Pharmacol.*, 65: 1701-1725. (2013)

TP-17 / SIMULTANEOUS DETERMINATION OF ARIPIRAZOLE, DEHYDRO-ARIPIRAZOLE, OLANZAPINE, RISPERIDONE, PALIPERIDONE, QUETIAPINE AND CLOZAPINE IN HUMAN PLASMA BY LC-MS/MS

Dóra Koller - Francisco Abad-Santos - Aneta Wojnicz

Hospital Universitario de la Princesa, Clinical Pharmacology Department, Universidad Autónoma de Madrid, Instituto de Investigación Sanitaria la Princesa, Madrid

Keywords: antipsychotics, LC-MS/MS, SPE, pharmacokinetics, quantification, method validation

Introduction

A simple and sensitive LC-MS/MS method was developed and validated for the simultaneous quantification of aripiprazole (ARI) and its active metabolite dehydro-aripiprazole (DARI), olanzapine (OLA), risperidone (RIS) and its active metabolite paliperidone (PAL), quetiapine (QUE) and clozapine (CLO) in human plasma. These drugs have been widely used in the treatment of schizophrenia and schizoaffective disorders [1-3].

Methods

Stable isotopically labeled internal standards - [$^2\text{H}_8$]-ARI, [$^{13}\text{C}, ^2\text{H}_3$]-OLA, [$^2\text{H}_4$]-RIS, [$^2\text{H}_4$]-PAL, [$^{13}\text{C}_4$]-QUE and [$^{13}\text{C}, ^2\text{H}_3$]-CLO - were used for all analytes. Three-step phospholipids-eliminating microelution-solid-phase extraction was used for analyte extraction. The chromatographic separation was performed during 8 min under gradient conditions through a mobile phase consisting of 0.2% formic acid in water, pH=3 and acetonitrile with a flow rate of 0.6 ml/min.

Results

The analytical method has been validated according to the recommendations of European Medicines Agency (EMA) and Food and Drug Administration (FDA) through tests of precision, accuracy, stability, sensitivity and specificity. The coefficient of variation was not more than 20% for the lower limit of quantification and was below 15% for the rest of quality controls. The method is able to detect the drug levels in the ranges of 0.18-120 ng/ml for ARI, 0.25-80 ng/ml for DARI, 0.5-1000 ng/ml for CLO, 1-100 ng/ml for OLA, 0.7-60 ng/ml for RIS, 0.2-30 ng/ml for PAL and 0.5-160 ng/ml for QUE during therapeutic drug monitoring. Comparing with protein precipitation (PPT), our SPE sample preparation method resulted to be

effective in removing phospholipids from plasma.

Conclusions

We have developed a sensitive LC-MS/MS method to determine the levels of 5 antipsychotics and 2 active metabolites with the usage of stable isotopically labeled internal standards in human plasma. The method is currently being utilized in pharmacokinetic studies. Its application could help to optimize the effectiveness of treatment and manage side effects or toxicity leading to individualized therapeutic decisions for each patient.

Novel Aspect

The present method applies effective phospholipids' removal during SPE and therefore results in exceptional recoveries and low matrix effect.

References

1. C. Wang, W. Shi, C. Huang, J. Zhu, W. Huang, y G. Chen; Ann. Gen. Psychiatry, 16, 47, 2017.
2. D. Pilon et al. ; J. Med. Econ., 21, 135-143, 2018.
3. F. Franza, B. Solomita, G. Aldi, y G. Del Buono; Psychiatr. Danub., 29, n.o Suppl 3, 405-408, 2017.

TP-18 / LC-MS IDENTIFICATION AND QUANTIFICATION OF BLEOMYCININ IN SERUM AND TUMOR TISSUE

Tina Kosjek⁽¹⁾ - ***Tjasa Gornik***⁽¹⁾ - ***Dusan Zigon***⁽¹⁾ - ***Masa Bosnjak***⁽²⁾ - ***Simona Kranjc***⁽²⁾ - ***Mojca Krzan***⁽³⁾ - ***Ales Groselj***⁽⁴⁾ - ***Maja Cemazar***⁽²⁾ - ***Gregor Sersa***⁽²⁾

Jozef Stefan Institute, Dept. Envir. Sci., Ljubljana (1) - Institute of Oncology, Department of Experimental Oncology, Ljubljana (2) - University of Ljubljana, Faculty of Medicine, Ljubljana (3) - University Medical Centre Ljubljana, Department of Otorhinolaryngology and Cervicofacial Surgery, Ljubljana (4)

Keywords: Bleomycin, metal complex, serum, tumor, LC-MS analysis

Introduction

Bleomycin (BLM) is a cytotoxic antibiotic that appears as a mixture of more than ten structurally strongly related glycopeptides, most abundant being BLM-A2 and BLM-B2 fractions. In vivo BLM chelates with several metals, in particular with copper. Aiming to investigate its biodistribution and pharmacokinetics as a part of a wider electrochemotherapy study, we developed a specific and sensitive mass spectrometry-based analytical method for determination of BLM in serum and tumor tissues.

Methods

Frozen tumors were grinded to fine powder, the suspension was sonicated, the supernatant filtered and loaded onto a preconditioned solid phase extraction (SPE) cartridge. Reversed-phase SPE of diluted serum and tumor suspensions was performed as the enrichment step. BLM behavior was investigated with UHPLC coupled to a QTOF MS[1]. To improve quantification characteristics the instrumental method was later translated to UHPLC coupled to QTRAP MS.

Results

The optimized sample preparation method used OasisHLB™ cartridges at acidic pH of the sample matrix, water as the wash solvent, and 40% methanol followed by acetonitrile for BLM elution. For UHPLC separation we selected BEH Amide stationary phase, which clearly and distinctively separated the Cu complexes of BLM-A2 and BLM-B2 fractions.

The investigation of BLM chelating properties implies that Cu complexes form in vitro and in vivo, and offer a viable form for the quantitative determination of BLM in biological samples. The amount of Cu available in the serum is sufficient to quantitatively chelate metal-free BLM even at its highest concentrations.

Metal-free BLM-A2 shows a 2⁺ charged fragment ion at *m/z* 707.76, but it is a "transient" species that forms

the Cu complex. The latter shows an abundant ion at m/z 707.21, and m/z 707.71 as the second isotope ion in the typical Cu isotopic pattern. In the complex mixture of several BLM fractions and their metal complexes, and considering also background ions from the biological tissue matrices, this second isotopic ion may easily be misinterpreted and incorrectly assigned to the metal-free BLM-A2, which could be the reason, why the attempts to determine BLM with MS have been so scarce.

Conclusions

Except in the controlled metal-free environment BLM is not found in a free form, but is highly prone to form chelates. Thus, MS analysis results in the apparent disappearance of BLM characteristic fragment ions. In addition, the chelates may be hidden at first sight, since they tend to form double charged ions. Choosing a HRMS proved crucial for recognizing BLM MS behavior, since it enabled us to follow the isotopic pattern of the double charged ions in the clusters.

Novel Aspect

The LC-MS analytical method for determination of bleomycin in serum and tumor tissues is developed for the first time.

References

1. Kosjek T., Krajnc A., Gornik T., Zigon D., Groselj D., Sersa G., Cemazar M., *Talanta* 160, 164–171 (2016).

TP-19 / SOLID PHASE MICROEXTRACTION FOR THE QUANTIFICATION OF DOXORUBICIN AND UNTARGETED METABOLIC PROFILING OF LUNG TISSUE DURING IN VIVO LUNG PERFUSION

Nikita Looby⁽¹⁾ - *Barbara Bojko*⁽²⁾ - *Germán Augusto Gómez-Ríos*⁽¹⁾ - *Marcos Tascon*⁽¹⁾ - *Anna Roszkowska*⁽³⁾ - *Krzysztof Goryński*⁽²⁾ - *Janusz Pawliszyn*⁽¹⁾

University of Waterloo, Department of Chemistry, Waterloo (1) - Collegium Medicum, Nicolaus Copernicus University in Torun, Department of Pharmacodynamics and Molecular Pharmacology, Bydgoszcz (2) - Medical University of Gdansk, Medical University of Gdansk, Department of Pharmaceutical Chemistry, Gdansk (3)

Keywords: Chemotherapy, in vivo lung perfusion, human trial, drug monitoring, solid phase microextraction.

Introduction

In vivo lung perfusion is a novel technique designed to localize chemotherapy to the lung tissue suffering metastasis. Effectively, increased doses of chemo agents can be used while reducing systemic exposure to the individual thus decreasing associated side effects (1). As there are currently no techniques available to determine drug tissue concentration during IVLP without taking biopsies, SPME emerged as a feasible strategy for in vivo quantification of doxorubicin (2).

Methods

In vivo lung perfusion (IVLP) was used to administer doxorubicin for treatment of metastasized lung tumor in a human being. Bio SPME fibers were inserted in vivo into three sections of the lung to monitor drug concentration, distribution and effect in lung tissue during IVLP. Liquid chromatography tandem mass spectrometry was used for drug quantification, while high resolution mass spectrometry was employed for untargeted analysis to determine markers of tissue stress or damage.

Results

An administered dose of 5 $\mu\text{g}/\text{mL}$ of doxorubicin in perfusate solution resulted in a drug tissue concentration 40 times higher than the administered dose, reaching a concentration of 200 $\mu\text{g}/\text{mL}$ within the first hour of IVLP. This increase was marked by a $\sim 50\%$ decrease in the perfusate drug concentration.

Subsequent analysis of the perfusate solution as well as the lung tissue at hourly intervals showed a decreasing trend in the concentration in both biological matrices. Untargeted metabolic profiling processed via multivariate analysis using Metaboanalyst 4.0 showed an interesting feature at m/z 166.0864 tentatively identified as L-phenylalanine, an endogenous metabolite that could serve as a possible marker of tissue stress or damage since it is known to be a metabotoxin at high concentrations.

Conclusions

Targeted analysis for doxorubicin quantitation and metabolic profiling to assess chemotherapeutic effect in lung tissue could be achieved using SPME. A decreasing trend of doxorubicin concentration in perfusate is linked to its continuous accumulation in tissue during IVLP. Furthermore, it is hypothesized that the decreasing doxorubicin concentrations in lung tissue may be due to intracellular uptake of the drug, making it less available to the SPME fibers that extract via free concentration.

Novel Aspect

For the first time, concurrent semi-quantification of doxorubicin and chemical tissue biopsies in vivo in human lung during surgery can be achieved via SPME.

References

- 1 P. Reck dos Santos, J. Sakamoto, M. Chen, V. Linacre, C. Arce, M. Liu, T. K. Waddell, S. Keshavjee, M. Cypel, *Ann. Thorac. Surg.*, 2016, 101, 2132–2140.
- 2 A. Roszkowska, M. Tascon, B. Bojko, K. Goryński, P. R. dos Santos, M. Cypel, J. Pawliszyn, *Talanta*, 2018, 183, 304–310.

TP-20 / ASSESSMENT OF ANTIBODY-DERIVED THERAPEUTICS AT THE INTACT AND MIDDLE-UP LEVEL BY CESI-MS

Elena Domínguez-Vega⁽¹⁾ - *Rob Haselberg*⁽²⁾ - *Govert W. Somsen*⁽²⁾ - *Jim Thorn*⁽³⁾ - *Christopher Lößner*⁽⁴⁾ - *Stephen Lock*⁽³⁾

Leiden University Medical Center, Center for Proteomics and Metabolomics, Leiden (1) - Vrije Universiteit, Division of BioAnalytical Chemistry, Amsterdam (2) - SCIEX, CE and Biopharma, Warrington (3) - SCIEX, CE and Biopharma, Darmstadt (4)

Keywords: CE-MS, CESI, antibody-based pharmaceuticals, nanobodies, proteoforms

Introduction

Antibody-based pharmaceuticals often contain complex structural heterogeneity which requires enhanced analytical methods for reliable characterization of variants and degradation products. In this poster we will describe how CE-MS (CESI-MS which integrates CE and electrospray ionization (ESI) into a single device) in combination with high-resolution MS detection has been used for profiling antibody therapeutics.

Methods

By using a neutral capillary coating (to provide near-zero electroosmotic flow) and an acidic background electrolyte intact model proteins were detected with overall migration-time RSDs below 2.2% (using 3 different capillaries). Samples, including mono- and bivalent nanobodies, and three monoclonal antibodies (mAbs) were tested.

Results

Intact nanobodies were resolved from their degradation products, which could be assigned to deamidated, cleaved, and truncated forms at the C-terminal tag with excellent resolution of isomeric deamidated products obtained.

The mAbs were analysed intact and after digestion by the IdeS endoproteinase (middle-up approach). CE-MS of intact mAbs resolved clipped species (*e.g.* light chain and light chain-heavy chain fragments) from the native protein as well as glycoforms containing sialic acids from their non-sialylated counterparts. For IdeS-digested mAbs, F(ab)₂ and Fc/2 were efficiently resolved. While migration of Fc/2 fragments were fairly similar for the three mAbs, the migration of the F(ab)₂ strongly depended on the mAb. All Fc/2 charged variants, which included glycoforms containing sialic acids and other PTMs such as loss of C-terminal lysine or deamidation of Asn, were nicely separated in less than 20 min. This allowed a detailed and reliable assessment of the Fc/2 heterogeneity (18-33 proteoforms) for the mAbs studied.

TP-21 / DETECTION AND CHARACTERIZATION OF STABLE IMINIUM ION REACTIVE METABOLITE USING HIGH RESOLUTION MASS SPECTROMETRY

Tommaso Miraval

Idorsia LTD, DMPK, Allschwil

Keywords: Drug Metabolism, Reactive Metabolites, High Resolution Mass Spectrometry

Introduction

Iminium ion reactive metabolites.

Xenobiotics including drugs may be metabolized to form reactive metabolites that covalently bind to protein or DNA. Such bioactivation is suspected to play a causative role in various adverse drug reactions. Thus, early detection of reactive metabolites has become an important part of the drug discovery process. Several tertiary amine-bearing drugs undergo P450- or MAO-catalyzed oxidation to unstable iminium ion species which readily hydrolyze to aldehydes, both of which being potentially reactive. However, due to their inherent chemical instability, the use of trapping agents (*e.g.* cyanide, methoxyamine) is required to put them into evidence. In this poster, we report unusual iminium ion metabolites which appear sufficiently stable to be directly detected and characterized using liquid chromatography (LC) coupled with high resolution mass spectrometry (HRMS).

Methods

In vitro incubations and metabolites identification.

10 µM incubations of compound P, bearing a piperazine moiety, were performed with human liver microsomes in potassium buffer (0.1 M, pH 7.4) for 60 minutes with and without 1 mM solution of potassium cyanide. The reaction was initiated by addition of NADPH after 5 minutes of preincubation and stopped with an equal volume of acetonitrile.

10 µM incubations of compound P were also performed using fresh rat and cryopreserved human hepatocytes in suspension at cell density of 1 million cell/mL. An aliquot of each sample was taken at various time points and reaction was stopped with an equal volume of acetonitrile. After centrifugation at 13000 RPM for 10 minutes, all the samples were transferred to an HPLC vial and submitted to the Q-Exactive HF (ThermoFisher AG, Reinach) LC-MS system. After a preliminary full scan experiment in positive ion mode, representative samples, using a blank matrix without compound as a control, were chosen and then processed using Compound Discoverer 2.0 (ThermoFisher AG, Reinach) to identify the main metabolites. After the selection of molecular ions of interest, production reaction monitoring (PRM) experiments were performed to elaborate the fragmentation pattern of each metabolite. To refine structural elucidation, a deuterium exchange experiment was also performed, in which the number of exchangeable proton present in the molecule were established.

Results

Metabolic Profiling.

Metabolic profiles were similar between human microsomes and hepatocytes. No conjugated phase II

metabolite was observed whereas parent drug (P) was the most prominent entity detected in either matrices.

Several prominent phase I metabolites were also seen. This included the *N*-dealkylated product M1, which was further oxidized to M2 and M3, the *N*-oxide M4 and a doubly dehydrogenated product M5. In addition, the dehydrogenated M7 and amide derivatives of M1 (M8 and M9) and of P (M10) were also detected in hepatocytes.

In the rat cells, the nitrene derivative of M1 was the main entity detected. Parent drug (P) M1, M2 and M3, M4, M5, M8, M9, M10 were also observed. Based on fragmentation patterns, all biotransformations were assigned to the piperazine ring.

Both M5 and M7 did not exhibit any shift in the deuterium exchange experiment, suggesting that these metabolites were already charged and did not require an ionization proton for MS detection. Metabolite M7 was successfully trapped with the iminium trapping agent cyanide in human liver microsomes, and the concomitant lowering of M5 abundance was observed.

Conclusions

The involvement of an intermediary iminium ion metabolite in the metabolic pathway of P was initially suspected from characteristic downstream metabolites of such species, the dealkylated M1 and the amide derivative M10 only seen in hepatocytes, which would be consistent with a common aldehyde oxidase-mediated detoxification process. This was further substantiated by the successful trapping of the postulated iminium ion in the presence of cyanide. However, the detection and characterization of M7, which appears to be the iminium ion itself, was quite unexpected. The potential conversion to a dimer or an enamine derivative were excluded based on MS and deuterium exchange data, respectively. The exact position of the unsaturated bond on the piperazine ring could not be allocated with more precision. The presence of the pyrazinium metabolite M10 suggested that M7 was also undergoing secondary metabolism in the incubation matrices.

Novel Aspect

Based on these results M7 was assigned as an unusually stable iminium ion species. This metabolite appeared not only stable in solution and during LC-MS analysis, but it could still be detected after several days of sample storage in autosampler. The isolation and further characterization of this iminium ion metabolite by alternative analytical techniques (*e.g.* NMR) would be required for definitive structural elucidation.

References

1. PeterlinMasic L., Current Drug Metabolism, 12, 35-50 (2011).
2. Kalgutklar A. S., Chemical Research in Toxicology, 30, 220-238 (2017).
3. Argoti D., Liang L., Conteh A. et al., Chemical Research in Toxicology, 18, 1537-1544

TP-22 / PHYTOCHEMICAL COMPOSITION OF EXTRACTS FROM WOOD AND LEAVES OF LEBANESE CEDAR (CEDRUS LIBANI) USING GC/MS AND LC/MS-MS

Samia Mokh⁽¹⁾ - **Mohamad Al Iskandarani**⁽²⁾ - **Marcelo Locatelli**⁽³⁾ - **Simone Carradori**⁽³⁾ - **Christian Celia**⁽³⁾ - **Jaber Farouk**⁽⁴⁾

Lebanese Atomic Energy Commission, Laboratory of Analysis of Organic Compound, Beirut, Lebanon (1) - Lebanese University, Faculty of Public Health, Laboratory, Beirut, Lebanon (2) - Università degli Studi G. d'Annunzio Chieti e Pesc., Department of Pharmacy, Chieti, Italy (3) - Laboratoire d'Analyse de Pesticides et de Polluants Organiques Commission Libanaise de l'Energie Atomique - CNRSL (4)

Keywords: Lebanese cedar, extraction methods, Characterization, GC/MS

Introduction

Lebanon Cedar is a native Mediterranean species with a pleasant odor and valuable wood. Despite the efforts devoted to the study of this plant, it remains little characterized from the analytical, technological and biological point of view. The aim of this study is to determine the phytochemical profile of leaves and wood of *cedrus libani* by using GC/MS and extract all components by maceration, soxhlet and ultrasound-assisted extraction (UAE) techniques.

Methods

Lebanese cedar leaves and wood samples were collected from a natural reserve in Nabatieh, southern Lebanon. These samples were extracted using three different methods: ultrasound assisted extraction, maceration and Soxhlet extraction. Four different extraction solvents were used: hexane, ethyl acetate, acetonitrile and methanol. The extracts thus obtained were analyzed by gas chromatography coupled with mass spectrometry (GC/MS).

Results

The results show that among the chemical classes extracted, five are interesting: terpenes, polyphenols, phytosterols, vitamins and cortisol. Hexane appears to be an ideal solvent for the extraction of terpenes. Methanol is an excellent solvent for the recovery of polyphenols. About the phytosterols and vitamins, the four solvents are able to extract them. As for cortisol, only methanol makes it possible to obtain it. On the other hand, ultrasound assisted extraction seems the best method of extracting the five chemical classes of interest. However, Soxhlet extraction due to the high level of hydrocarbons extract does not have the ability to use for the extraction of a very large number of compounds as compared to the other two methods. For maceration, its weak point is that it cannot extract the cortisol. In addition, several compounds of medical and industrial importance were extracted from lebanese cedar, such as vitamin E, vitamin A1, 24,25-dihydroxyvitamin D3, and phytol.

Conclusions

Analysis of these results shows that the ultrasound assisted extraction using methanol solvent of leaves and wood was the best method in order to obtain the five interesting families: terpenes, polyphenols, phytosterols, vitamins and cortisol. The extraction method was optimized and the identification of compound was completed by the GC/MS using the library of NIST and should be completed by LC/HRMS and LC/RMN by putting in evidence its analytical fingerprint.

Novel Aspect

The results illustrate the importance of optimizing the extraction of the principal families present in the leaves and wood of *libani cedrus* prior to their measurement using GC/MS.

References

1. Al-Mansoub, M., Asmawi, M., Murugaiya H., Journal of the Science of Food and Agriculture, 94(8), 1552–1558 (2014).
2. Asbahani, A. El, Miladi, K., Badri, W., Sala, M., Addi, E. H. A., Casabianca, H., Elaissari, A., International Journal of Pharmaceutics, 483, 220–243 (2015).
3. Azwanida N. N., Medicinal & Aromatic Plants, 4, 3–8 (2015).
4. Fleisher, A. & Fleisher, Z., Journal of Essential Oil Research, 12, 763–65 (2011).
5. Espada-Bellido, E., Ferreiro-González, M., Carrera, C., Palma, M., Barroso, C. G., Barbero, G. F., Food Chemistry, 219, 23–32 (2017).

TP-23 / A COMBINATION OF CHEMICAL PROTEOMICS AND DARST-BASED APPROACHES TO DISCLOSE NATURAL PRODUCTS BIOLOGICAL TARGETS

Federica Del Gaudio ⁽¹⁾ - **Sara Ceccacci** ⁽¹⁾ - **Federica Pollastro** ⁽²⁾ - **Matteo Mozzicafreddo** ⁽¹⁾ - **Alberto Minassi** ⁽²⁾ - **Raffaele Riccio** ⁽¹⁾ - **Maria Chiara Monti** ⁽¹⁾

Dept. of Pharmacy, University of Salerno, Fisciano (Salerno) (1) - Dept. Pharmaceutical Sciences, Università degli Studi del Piemonte Orientale Amedeo Avogadro, Novara (2)

Keywords: target identification, affinity based chemical proteomics, DARTS (drug affinity responsive target stability), natural molecules.

Introduction

Disclosing Natural Products Interactome.

The identification of the biological and potential therapeutic targets of small molecules, along with their off-targets, is a significant concern. Here, a combination of the affinity-based chemical proteomics (AP-MS) and a simple, universally applicable approach based on the drug affinity responsive target stability (DARTS) has been applied to disclose the interactome of Arzanol, a natural prenylated phloroglucinylypyrone, used in folk medicine to treat infections and inflammatory diseases [1].

Methods

Combining AP-MS and DARST approaches

AP-MS [2] is based on the preparation of an Arzanol bearing solid support, the isolation of its potential targets from a crude cell extract, SDS-PAGE of Arzanol bound targets and their identification by MS. The DARTS approach [3] takes advantage of a reduction in the protease susceptibility of Arzanol targets, easily revealed by SDS-PAGE-LC-MS analysis or by immunoblotting, without Arzanol modification or immobilization.

Results

Brain Glycogen Phosphorylase as a novel Arzanol target.

AP-MS revealed brain glycogen phosphorylase (bGP) as the major target of Arzanol together with the muscle GP isoform and ATP-dependent 6-phosphofructokinases. DARTS experiments followed by both nano-LC-MS and immunoblotting confirmed bGP as a direct interactor of Arzanol. bGP catalyzes the rate-limiting step of glycogen mobilization in brain; its allosteric modulation by AMP, ADP, ATP and phosphorylation, finely tunes its activity in response to different signals. Thanks to AP-MS competition assays, the interaction site of Arzanol on bGP has been disclosed.

Conclusions

Expanding Arzanol potential bioactivity

The natural heterodimeric polyphenol Arzanol is known for its anti-inflammatory and antioxidant properties [4]. Up to now, a full screening of its biological most relevant interactors was still missing and, thus, a combination of chemical proteomics and DARTS experiments has been carried out on this intriguing natural compound, disclosing its interaction with bGP. Brain glycogen metabolism plays a critical role in major brain functions such as learning or memory consolidation: alteration of glycogen metabolism and glycogen accumulation in the brain contributes to neurodegeneration

Novel Aspect

DARST as the ideal approach to support AP-MS data.

In the past, AP-MS has been widely employed to identify targets of chemically reactive drugs. DARTS strategy is an ideal alternative and/or complement to AP-MS working on native unmodified molecules.

References

1. Appendino G., Ottino M., Marquez N., Bianchi F., Giana A., Ballero M., Sterner O., Fiebich B. L., Munoz E. *J Nat Prod.*, 70, 608 (2007)
2. Cassiano C., Esposito R., Tosco A., Zampella A., D'Auria M. V., Riccio R., Casapullo A., Monti M. C. *Chem Commun.*, 50, 406 (2014)
3. Lomenick B., Hao R., Jonai N., Chin R. M., Aghajan M., Warburton S., Wang J., Wu R. P., Gomez F., Loo J. A., Wohlschlegel J. A., Vondriska T. M., Pelletier J., Herschman H. R., Clardy J., Clarke C. F., Huang J., *PNAS*, 106, 21984 (2009)

4. Rosa A., Deiana M., Atzeri A., Corona G., Incani A., Melis M. P., Appendino G., Dessì M. A. *Chem Biol Interact.* 165, 117 (2007)

TP-24 / FOCAL ADHESION KINASE BINDING PARTNER STUDY IN HUMAN COLON CANCER CELL HCT-116 USING IMMUNOPRECIPITATION METHOD COUPLED WITH NANO-LC/MSMS TECHNIQUE

Binh Thanh Nguyen - Min-Jung Kang

Korea Institute of Science and Technology, Korea University of Science and Technology, Seoul

Keywords: Binding partners, Colon cancer cell HCT-116, Focal Adhesion Kinase, Immunoprecipitation, Proteomics

Introduction

Focal adhesion kinase (FAK) is a non-receptor tyrosine that regulates migration, proliferation and survival of cells. A large number of studies have strongly implicated FAK in the development of colon cancer. The novel potential functions of FAK in the regulation of epithelial-mesenchymal transition (EMT) has been studied extensively in invasive and metastatic cancers. The study focuses on new binding partners of FAK as therapeutic targets for colon cancer.

Methods

Human colon cancer HCT-116 was cultured, harvested and lysed. The protein concentration was measured by BCA kit. FAK expression was confirmed by Western blot. FAK and binding partners were captured by immunoprecipitation method and separated by SDS-PAGE. Gels were excised and digested with trypsin. Tryptic peptides were injected to nano-LC/MSMS system and data were searched against the Mascot search engine. The binding candidates were analyzed with Western blot.

Results

FAK was expressed highly in colon cancer cell HCT-116 compared to the other cancers and normal colon cell. Immunoprecipitation of FAK protein was optimized in the manner of bead selection, cell lysate amount, antibody amount, washing conditions, and elution volume. Gel-LC/MSMS data and mascot search revealed about 350 protein identity. FAK was found with the highest score indicating the successful pull-down and digestion. Among the matched proteins, reportedly known binding partners of FAK include growth factor receptor-bound protein 14 GRB14, integrin-linked protein kinase ILK, Ras-related protein RAB-22A, Rho-related GTP-binding protein RhoU. New candidates for FAK binding partners proposed in our result include SACS, saccin; SYNE3, nesprin; ZYX, zyxin; PLK1, serine/threonine protein. The new binding partner candidates were confirmed using Western blot.

Conclusions

The further work include the effect of growth factors, FAK inhibitors on colon cancer cells to elucidate the upstream regulators or downstream regulators involving the FAK signaling pathways. New binding partners of FAK once validated and characterized will be potential biomarkers for colon cancer diagnosis and useful for the development of therapeutic targets.

Novel Aspect

The literatures have not shown evidences for interaction between new binding candidates of FAK. Necessary confirmation would be novel for therapeutic target development.

References

1. K.M. Quesnelle, A.M. Sparano, Y. Wang, G.S. Weinstein, M.S. Brose, *Cancer research*, 66, 798-798 (2006)
2. J. E. Hall, W. Fu, M. D. Schaller, *International review of cell and molecular biology*, 288 185-225(2011)
3. S. T. Arold, M. K. Hoellerer, M. E. M. Noble, *Structure*, 10 319-327 (2002)

4. J. Zhao, J. L. Guan, *Cancer Metastasis Reviews*, 28 35-49 (2009).

TP-25 / COMPREHENSIVE N-GLYCAN STRUCTURE CATALOGUE OF THERAPEUTIC ANTIBODIES

Myung Jin Oh - Nari Seo - Hyun Joo An

AGRS, Chungnam National University, Daejeon

Keywords: monoclonal antibody, *N*-glycan, reference, mass spectrometry

Introduction

Glycosylation of monoclonal antibody (mAb) therapeutics plays vital roles in drug's biological activity and safety. Monitoring of glycosylation affected by various environment conditions is highly required to access QA/QC of mAbs. Here, we developed a LC/MS based *N*-glycan library for rapid and high-throughput identification in various mAbs.

Methods

The *N*-glycan catalogue was constructed using human immunoglobulin G and five different mAb therapeutics produced by various cell-based expression systems. Reduced native *N*-glycans were separated and analyzed by porous graphitized carbon (PGC)-chip nano-LC/MS. Accurate masses were used to compositionally profile *N*-glycans and their structures were further elucidated by CID MS/MS. Major glycans (G0F, G1F, and G2F) on mAbs were quantified by UHPLC-HILIC-FLD.

Results

Representative commercial mAbs (adalimumab, bevacizumab, infliximab, rituximab, and trastuzumab) produced by CHO cell and Sp2/O cell were selected to construct detailed *N*-glycan catalogue for rapid characterization. In parallel, human IgG was used as the reference to confirm human-like glycan structures. The majority of the *N*-glycans from five commercial mAbs were characterized by matching the RT and accurate masses with the structure of glycans found in a human IgG reference. Corresponding compositions were further confirmed using structural connectivity of the individual monosaccharide residues by tandem MS. In particular, non-human glycans found in non-human mammalian cell system were newly identified by accurate masses, RTs, and tandem MS. Antibody *N*-glycan catalogue containing over 70 distinct glycan compositions provides structures, relative quantity, and chromatographic retention time (RT) based on PGC column.

Conclusions

We created the *N*-glycan catalogue for rapid identification of glycans in mAbs. The combination of accurate masses and LC retention times enables real-time identification for glycan structures. Product ions of each *N*-glycan by tandem MS were also recorded in the catalogue. *N*-glycan catalogue of mAbs can be used as a valuable reference in order to evaluate mAb variants and/or biosimilars for both developmental and regulatory purposes.

Novel Aspect

Comprehensive *N*-glycan catalogue makes it possible to identify their structures in a high-throughput manner by matching LC retention times and accurate masses.

TP-26 / MONITORING OF SEQUENCE VARIANTS BY MAM USING HIGH RESOLUTION MASS SPECTROMETRY

***Kerstin Pohl*⁽¹⁾ - *Yihan Li*⁽²⁾ - *Chris Nortcliffe*⁽³⁾ - *Annu Uppal*⁽⁴⁾ - *Ferran Sanchez*⁽⁵⁾ - *Jörg Dojahn*⁽¹⁾ - *Sibylle Heidelberger*⁽³⁾**

SCIEX, Biologics, Darmstadt (1) - SCIEX, Biologics, Redwood City (2) - SCIEX, Biologics, Warrington (3) - SCIEX, Biologics, Gurugram (4) - SCIEX, Biologics, Barcelona (5)

Keywords: Biologics, Sequence Variants, Data-Independent, Antibodies, MAM

Introduction

Sequence variants (SVs) are unintentional amino acid substitutions and contribute to the overall heterogeneity of biotherapeutics. The presence of these replacements can be at differing levels, dependent on the cause of the variant, making identification difficult. Here, we monitor the low level SVs using a highly sensitive, high resolution mass spectrometer, the SCIEX TripleTOF® 6600, coupled to built-for-purpose BioPharmaView™.

Methods

NIST mAb reference standard (#RM8671) was digested using standard conditions. First, digested NIST was run using a 3 hour gradient on a ExionLC™ AD system coupled to a TripleTOF® 6600 and data was acquired in DDA mode. Data was searched using ProteinPilot™ including a search for SVs. Subsequently NIST was acquired using SWATH® acquisition. The data was processed in BioPharmaView™ software to identify and relatively quantify the sequence variants.

Results

NIST mAb was digested and run using a 3 hour reversed-phase gradient in data-dependent acquisition mode. The data was searched with ProteinPilot™ with an emphasis on IgG antibody and SV identification in order to look for variants. Once the SVs were identified, these sequences were input into BioPharmaView™ as impurities to monitor with SWATH® acquisition on a shorter gradient. Using the peptide mapping workflow including MAM analysis in BioPharmaView™, SVs can be identified and also quantified based on the non-variant peptide despite the shorter gradient and their low abundance. The data-independent SWATH® acquisition ensures that MS/MS data is available on the whole mass range throughout the gradient.

SVA requires a mass spectrometer with a high linear dynamic range, high sensitivity and speed to ensure that the variants are seen at MS level and also fragmented resulting in a good quality MSMS data.

Conclusions

SVs contribute to the overall heterogeneity of biologics and are a measure of product, and production stability along with potential detrimental effects on physicochemical properties, tertiary structure, and efficacy of the product. Misfolded proteoform variants may be toxic impurities. The ability to detect these variants at low level assists manufacturers in determining the potential consequences of the variant and differences from the expected product.

Novel Aspect

Sequence variant analysis using MAM workflow.

TP-27 / EPI TOPE PEPTIDES IDENTIFIED BY ONLINE BIOSENSOR- MS NEUTRALIZE PATHOPHYSIOLOGICAL ANTIBODIES AND OPEN NEW CLINICAL THERAPY APPROACHES FOR LYSOSOMAL DISEASES

Michael Przybylski⁽¹⁾ - *Loredana Lupu*⁽¹⁾ - *Pascal Wiegand*⁽¹⁾ - *Zdeněk Kukačka*⁽¹⁾ - *Stefan Maeser*⁽¹⁾ - *Yannick Baschung*⁽¹⁾ - *Friedemann Völklein*⁽²⁾ - *Alexander Lazarev*⁽³⁾ - *Anna Maria Papini*⁽⁴⁾ - *Julia Hennermann*⁽⁵⁾ - *Katja Schmitz*⁽⁶⁾

Centre for Biopolymer Analysis and Biomedical Mass Spectrometry, Rüsselsheim am Main (1) - Rhein Main University, Rüsselsheim am Main (3) - Pressure Biosciences Inc., Pressure Biosciences Inc., Boston (4) -

University of Florence, Florence (5) - Universitätsmedizin Mainz, Mainz (6) - Techn. University Darmstadt, , Darmstadt (7)

Keywords: Biosensor-MS combination; Antibody epitope identification; Lysosomal diseases; Enzyme replacement therapy; Neutralizing epitope peptides.

Introduction

Epitope analysis is a key step in the development of therapeutic antibodies. Lysosomal Storage Diseases (LSDs) are successfully treated by Enzyme Replacement Therapy (ERT); however substantial, life threatening effects arise from antibodies that block the infused enzyme. A new online proteolytic excision SPRMS combination (PROTEX-MS) enables simultaneous structure identification and affinity quantification of antibody epitopes. Here we report a new therapeutic concept by identification of antibody epitopes upon ERT. Synthetic epitopes against pathophysiological antibodies open new clinical therapy approaches by (i), hyposensitization of patients, and (ii) molecular apheresis to deplete antibodies.

Methods

Epitopes from antisera of Fabry disease (FD; α -Galactosidase deficiency) patients were identified using the PROTEX-SPR-MS combination. Tryptic peptide mixtures were loaded onto the SPR-MS; after washing out nonbinding peptides, epitope peptides were eluted with 0.01% TFA into the affinity interface, ESIMS was performed with a Waters Quattro Ultima MS. The epitope peptide, α -Gal(309-332) was identified from 2 FD patients and revealed high affinity (KD, 39 nM) comparable to the full length enzyme (KD, 16 nM). Linear, stabilized epitope peptides were prepared by SPPS.

Results

Key tools of the online PROTEX SPR-MS Epitope Analyzer are (i), a new SPR suitable for direct application of ESI-MS without any buffer change (T-SPR); (ii), a new high pressure proteolytic system with substantially enhanced efficiency of protease digestion. The PROTEX SPR-MS analyzer was applied to the epitope elucidation and affinity characterization of antibodies upon ERT in FD patients, and provided the epitope, α Gal(309-332), in which the lysine residues K-315 and 326 were shielded upon antibody binding. The epitope α Ga(309-332) was synthesized by solid phase peptide synthesis (SPPS) and purified by HPLC on a semipreparative C4-column, yielding homogeneous peptide with >95 % purity. SPR determination of the epitope peptide and the full length enzyme provided high affinities (KD, 39 and 16 nM, respectively). The α Gal(309-332) epitope contained two internal lysine residues that were not cleaved upon proteolytic epitope excision, consistent with the shielding upon antibody binding.

Conclusions

The clinical application potential of the PROTEX SPRMS Epitope Analyzer was shown in epitope elucidation from therapeutic and pathophysiological antibodies of lysosomal enzymes, with affinity constants (KD) amenable from milli- to nanomolar ranges. Epitope peptides in FD patients are currently developed for antibody neutralization in clinical ERT by systemic application and apheresis at physiological conditions. Thus, the identification, chemical synthesis, and biochemical evaluation of antibody epitopes open provide new clinical treatment approaches for reversing the immunogenicity and reconstituting the therapeutic efficacy of ERT in lysosomal diseases.

Novel Aspect

Antibody epitope identification in ERT of lysosomal diseases opens new clinical therapy to neutralize antibodies and reconstitute therapeutic efficiency.

References

1. Przybylski M., Lupu L., Voelklein F. (2018) Eur. Pat. Appl.
2. Iurascu M., Belaunzar O.M., Cozma C., Petrusch U., Renner C., Przybylski M. (2016) J. Am. Soc. Mass Spectrom. 27, 1105-1112.

3. Moise A., Maeser S., Rawer S., Eggers F., Murphy M., Bornheim J., Przybylski M. (2016) *J. Am. Soc. Mass Spectrom.* 27, 1071-1078.
4. Paraschiv G., Vincke C., Czaplewska P., Manea M., Muyldermans S., Przybylski M. (2013) *J. Mol. Recogn.* 26, 1-9
5. Kukacka Z., Iurascu M., Lupu L., Rusche H., Murphy M., Altamore L., Borri F., Maeser S., Papini A.M., Hennermann J., Przybylski M. (2018) *ChemMedChem.*, DOI 10.1002/201800094.

TP-28 / CHROMATOGRAPHY COUPLED TO MASS SPECTROMETRY FOR STABILITY STUDIES IN DRUG DISCOVERY

Laurence Queguiner - David Speybrouck - Sebastien Thomas

Janssen Research & Development, Discovery Sciences, Val de Reuil

Keywords: Mass Spectrometry, Chromatography, Stability studies, Pharmaceuticals

Introduction

Pharmaceutical companies always try to reduce the time spent on Drug Discovery to push forward compounds of interest into Clinical Development. To achieve the pre-clinical phase, Analysts need to provide the right and critical physicochemical data to the chemists.

Thanks to hyphenated techniques, we have developed in house-procedures and methods to test the stability of compounds, only looking for information which are relevant for the develop ability of a compound.

Methods

Depending on the compounds to be tested, we develop methods in SFC/MS for the chiral stability if compounds are pure enantiomers, after a screening to define the best separation. Chemical stability is monitored in UHPLC/MS either in LR or HRMS. If a compound is found unstable, High-Resolution equipment allows the identification of degradation compounds.

Our protocols are set up after discussion between chemists and analysts depending on the information they need.

Results

The chemical stabilities are monitored in UHPLC/UV/MS, on a UPLC/UV/MS (Quattro-micro) or a UHPLC/UV/MS (Q-Exactive) instruments, UV results are reported in tables and peak areas are compared to areas obtained at the beginning of the study, this is useful to discriminate precipitation and degradation. If degradation occurs, MS data are examined to give information about the degradation compounds to the chemists, if needed, HRMS data (full MS and after fragmentation) can provide structural information.

For the chiral stability, the same solutions are run in SFC/MS. SFC/UV having a low sensitivity, solutions are analyzed in MS in SIM (Selected Monitoring ions) mode. Most of the time, this technique allows to detect less than 1% of diastomer but to improve the sensitivity on poorly ionizable compounds or to discard salts coming from the highly charged buffers a step of Solid Phase Extraction (SPE) on cartridges can be added.

Conclusions

To push forward or discard compounds showing good or bad pharmaceutical potencies, the combination of different hyphenated techniques available in our laboratory helps the analytical chemists to fully describe the stability of a compound, depending on the product itself and the information needed by the chemists.

Novel Aspect

The combination of MS coupled to different chromatographic techniques are useful to obtain critical physicochemical data of pharmaceutical compounds.

TP-29 / ABCs OF mABs AND ADCs

Shaun Reeksting⁽¹⁾ - ***Terry Chapman***⁽²⁾ - ***Bayan Alkhawaja***⁽³⁾ - ***Andrew Watts***⁽³⁾ - ***Anneke Lubben***⁽¹⁾

University of Bath, Faculty of Science, Chemical Characterisation and Analysis Facility (CCAF), Bath (1) - Bath ASU, Research and Development laboratory, Corsham (2) - University of Bath, Pharmacy and Pharmacology, Bath (3)

Keywords: LC-MS, intact monoclonal antibody (mAB), antibody drug conjugate (ADC), drug-to-antibody ratio (DAR)

Introduction

Monoclonal antibodies (mABs) and antibody drug conjugates (ADCs) are a fast growing class of therapeutics. ADCs and mABs are complex heterogeneous populations of molecules. The NIST IgG mAB was used to validate our Agilent LC-nanoESI-CHIP-MS system performance. We additionally monitored the structural stability and degradation products of several intact and reduced forms of the Infliximab, Panitumumab, Rituximab, Kadcylla (ADC) and others.

Methods

For de-salting Amicon ultra 100 kDa spin filters, reduction with 4.4 molar eq. TCEP. For LC-MS Agilent CHIP cube on Agilent 6520, sample captured (0.2 µL inj on 40 nL Zorbax 300SB-C8, 5 µm) using 4µL/min H₂O + 0.1% formic acid (FA), then eluted off trap column onto the analytical column (43 mm x 75 µM with Zorbax 300SB-C8, 5 µm) using ACN:H₂O 90:10 with 0.1% FA. Processing with Agilent MassHunter MaxEnt deconvolution with mass using NIST Mass and Fragment Calculator, v1.32 [1].

Results

Evaluation of the system using the NIST IgG showed expected glycoforms for the intact protein [2]. Overall the mass accuracy for both intact and reduced mABs and ADC were within acceptable ranges (-12.1 and 20 ppm, respectively). For the Kadcylla ADC heterogeneous light chain D0, D1 and D2 proteins were shown with distinctive 957 Da emtansine payloads [3, 4]. In cases where TCEP was used to reduced mABs, partial or uncompleted reduction of disulfides was observed. This could have been due to potential steric hindrance of TCEP and inaccessibility of sequestered hydrophobic disulfides [5]. The inter-month mass stability and accurate mass reproducibility of an ADC suggested experiments should be run alongside controls for comparison. For ADCs under alkaline conditions and elevated temperatures the drug payload distribution was shown to be affected.

Conclusions

LC-nanoESI-MS is a useful analysis method, amongst other orthogonal techniques, to determine heterogeneous distributions of intact (non-reduced) and reduced forms of both mABs and ADCs.

Novel Aspect

LC-nanoESI-CHIP-MS system was suitable for intact and reduced analysis of mABs and ADCs and able to monitor drug payload stability on ADCs under accelerated stress conditions.

References

1. Eric L. Kilpatrick E., Liao W., Camara J., Turko I., Bunk D., Protein Expression and Purification, 85, 94-99 (2012)
2. Formolo T., Ly M., Davis D., Schiel J., State-of-the-Art and Emerging Technologies for Therapeutic Monoclonal Antibody Characterization Volume 2. Biopharmaceutical Characterization: The NISTmAb Case Study. 1-62 (2015)
3. Kadcylla / DRUG: Trastuzumabsequences, http://www.genome.jp/dbget-bin/www_bget?D03257+D09980

- Chen Y., Kim M., Zheng L., Deperalta G., Jacobson F., *Bioconjugate Chemistry*, 27, 2037–2047 (2016)
- Kirley T., Greis K., Norman A., *Biochemical and Biophysical Research Communications*, 480, 752–757 (2016)

TP-30 / DEVELOPMENT OF A HYBRID LC-MS/MS METHOD FOR THE QUANTIFICATION OF ANTIBODY-CONJUGATED PAYLOAD IN IN-VIVO AND IN-VITRO ADC SAMPLES

Federico Riccardi Sirtori⁽¹⁾ - *Elisa Bertotti*⁽¹⁾ - *Alessandra Ariaudo*⁽¹⁾ - *Francesco Molinaro*⁽¹⁾ - *Ilse De Salve*⁽¹⁾ - *Luca Barbero*⁽¹⁾ - *Seema Kumar*⁽²⁾ - *Simona Riva*⁽¹⁾ - *Theresa Goletz*⁽²⁾

RBM-Merck group, NBE Drug Disposition, QPD, Colletterto Giacosa (TO) (1) - EMD Serono, NBE Drug Disposition, QPD, Billerica, MA (2)

Keywords: Hybrid Method, LC-MS/MS, ADC, Antibody Conjugated Payload, Bioanalysis, NBE

Introduction

Antibody Drug Conjugates (ADC) are a novel class of targeted therapeutics that consist of a potent cytotoxic drug (payload) covalently bound to a monoclonal antibody (mAb) through a cleavable or non-cleavable linker [1,2]. ADCs can selectively deliver potent drugs through an antigen-antibody recognition mechanism. This step is followed by ADC-receptor complex internalization, enzymatic or chemical cleavage of the linker or protease degradation of Ab to intracellularly release the drug. The goal of ADCs is to improve therapeutic window by avoiding limitations related to classical antibody (Ab) based therapeutics or chemotherapy that uses potent cytotoxic small molecules.

Due to their multi-component structure and inherently heterogeneous nature, both antibody related analytes (total antibody, conjugated antibody) as well as small molecule related analytes (free payload, antibody-conjugated payload and total payload) are needed to understand ADC disposition in vivo. The bioanalytical assays to support ADC pharmacokinetics studies thus often require a combination of antibody (ligand binding assay, LBA) and small molecule analysis (LC-MS/MS) platforms.

Antibody-conjugated payload analyte describes the payload that is still conjugated to the antibody. Recently, there is a growing attention in Ab-conjugated payload determination to better define the ADC PK profiles [3,4,5].

Methods

A specific LC-MS/MS method was established and qualified for Ab-conjugated payload analyte for an ADC bearing a payload covalently bound through a cleavable linker. The method consists of an immunocapture of the ADC on a spin plate containing Protein A that selectively binds IgG antibodies. Subsequently a hydrolysis step was applied to the immobilized ADC on the plate to generate free payload. Eluates from the plate were subjected to a protein precipitation technique (PPT) step and supernatants are finally analyzed by LC-MS/MS to quantify the free toxin.

Results

The analytical method was qualified for linearity, accuracy, precision, selectivity, carryover, stability (freeze/thaw, benchtop, autosampler), and for dilution range.

To verify the reliability of the obtained results, a set of samples previously tested using conjugated payload method, were assayed by an orthogonal conjugated antibody assay, successfully qualified. This protocol, exploiting Meso Scale Discovery (MSD) technology, was developed and qualified to quantify ADC molecules in animal species (mouse, rat and cynomolgus plasma). A goat human IgG-heavy and light chain, monkey-adsorbed biotin labeled antibody was used as a capture antibody and a Sulfo-tag labeled monoclonal antibody against the toxin (was used as detection. The biotin labeled ADC target receptor was used as capture reagent instead of the goat anti human IgG to quantify the ADC in human plasma.

Conclusions

A good agreement between LBA conjugated Ab and LC-MS/MS conjugated payload results was found confirming the reliability of the two assays. The conjugated payload method was applied to support in-vitro (plasma stability studies) and in-vivo (PK) studies.

Novel Aspect

A novel hybrid LC-MS/MS method was developed and applied for the determination of Antibody Conjugated Payload in ADC bioanalysis. The procedure consists of an immunocapture step, digestion and detection of free toxin using an LC-MS/MS method.

References

1. Trail P. A., King H. D., Dubowchik G. M., *Cancer Immunology, Immunotherapy*, 52, 328–333 (2003).
2. Junutula J. R., Raab H., Clark S., Bhakta S., Leipold D. D., Weir S., *Nature Biotechnology*, 26, 925–932 (2008).
3. Sanderson R. J., Nicholas N. D., Lee C. B., Shawna Mae Hengel S. M., Lyon R. P., Benjamin D. R., C Alley S. C., *Bioanalysis*, 8, 55–63 (2016).
4. Liu A., *Journal of Chromatography B*, 1002, 54–62 (2015).
5. Xu L., Packer L. E., Qian M. G., Wu J. T., *Journal of Pharmaceutical and Biomedical Analysis*, 128, 226–235 (2016)

TP-31 / EPITOPE AND AFFINITY DETERMINATION OF RECOMBINANT MYCOBACTERIUM TUBERCULOSIS AG85B ANTIGEN TOWARDS ANTI-AG85 ANTIBODIES USING PROTEOLYTIC-AFFINITY MASS SPECTROMETRY AND BIOSENSOR ANALYSIS

Francesca Rinaldi⁽¹⁾ - **Loredana Lupu**⁽²⁾ - **Hendrik Rusche**⁽²⁾ - **Zdeněk Kukačka**⁽²⁾ - **Sara Tengattini**⁽¹⁾ - **Roberta Bernardini**⁽³⁾ - **Luciano Piubelli**⁽⁴⁾ - **Teodora Bavaro**⁽¹⁾ - **Stefan Maeser**⁽²⁾ - **Loredano Pollegioni**⁽⁴⁾ - **Enrica Calleri**⁽¹⁾ - **Michael Przybylski**⁽²⁾ - **Caterina Temporini**⁽¹⁾

University of Pavia, Department of Drug Sciences, Pavia (1) - Steinbeis Centre for Biopolymer Analysis and Biomedical Mass Spectrometry, Rüsselsheim am Main (2) - University of Rome "Tor Vergata", Department of Biology and Animal Technology Station, Rome (3) - University of Insubria, Department of Biotechnology and Life Sciences, Varese (4)

Keywords: epitope identification; affinity mass spectrometry; MALDI-TOF MS; SPR; glycoconjugate vaccine characterization

Introduction

Tuberculosis (TB) is the first cause of death from infectious diseases worldwide. Only a single anti-TB vaccine is currently available for clinical use, but its efficacy is not achieved with certainty [1]. The aim of this work is to provide a basis for the rational design of a neo-glycoconjugate vaccine against TB.

Methods

All proteins, glycoconjugates and antibodies used in this study were characterized by MALDI-TOF MS. The antibody epitope of recombinant Ag85B antigen was determined by proteolytic extraction-MS. Antibodies from sera of TB-patients, vaccinated subjects and a healthy control were employed, together with a commercial monoclonal antibody. Affinity binding constants were defined by SPR biosensor analysis.

Results

The MS epitope determination and SPR analysis allowed the identification and characterization of the interactions between Ag85B antigenic protein from Mycobacterium tuberculosis (MTB) and antibodies from different sources, both qualitatively and quantitatively.

It was possible to use human clinical samples to determine the Ag85B epitope. The identification of the

same binding areas of the protein interacting with a monoclonal antibody and with antibodies from human sera allowed to validate the simplified monoclonal model, confirming its applicability in SPR experiments. Therefore, an alternative screening method for the immunological evaluation of antigenic proteins from MTB has been developed. With this in vitro approach, the effects of the introduction of conservative mutations in Ag85B protein, as well as Ag85B conjugation with saccharidic moieties, were studied.

Conclusions

The combination of two analytical approaches yielded the identification of different assembled epitope regions on recombinant MTB antigens, their affinity binding constants in the interactions with specific antibodies and revealed the importance of protection from excessive glycosylation. These findings provide a molecular basis for an improved rational design of a new potential glycoconjugate vaccine against TB.

Novel Aspect

A qualitative and quantitative investigation of antigen-antibody interactions was performed to characterize a potential vaccine against TB.

References

1. World Health Organization, Global tuberculosis report 2017, Geneva: WHO Press (2017).

For information please contact: francesca.rinaldi01@ateneopv.it

TP-32 / MONITORING OF THERAPEUTIC DOSES OF ANTIPSYCHOTIC DRUGS IN SALIVA AND SERUM SAMPLES WITH THE USE OF SOLID-PHASE MICROEXTRACTION TECHNIQUE

Anna Roszkowska⁽¹⁾ - ***Lucyna Konieczna***⁽¹⁾ - ***Wioletta Żukowicz***⁽²⁾ - ***Marta Galewska***⁽¹⁾ - ***Wiesław Cubała***⁽³⁾ - ***Tomasz Bączek***⁽¹⁾

Department of Pharmaceutical Chemistry, Medical University of Gdansk, Gdańsk (1) - Department of Neurology, Medical University of Gdańsk, Gdańsk (2) - Department of Adult Psychiatry, Medical University of Gdańsk, Gdańsk (3)

Keywords: solid-phase microextraction, atypical antipsychotic drugs, therapeutic drug monitoring, LC-MS/MS

Introduction

The analysis of biological fluids to directly control the psychoactive substances is especially important in patients vulnerable to drug overdosing, e.g. pediatric and geriatric patients [1]. To provide fast and reliable information about the drug concentration, solid-phase microextraction (SPME) is introduced as a high-throughput sample preparation technique that facilitates accurate quantitation of the level of drugs in biological matrices [2,3].

Methods

Under optimized SPME conditions, the extraction of atypical antipsychotic drugs (risperidone, olanzapine, clozapine, quetiapine and aripiprazole) and internal standard (promethazine) from saliva and serum samples was performed with the use of C18 blades for 60 min at 1000 rpm. Then SPME blades were rinsed with nanopure water for 10 sec, and next desorption was performed for 10 min at 1000 rpm. Finally obtained extracts were analyzed with the use of LC-MS/MS system.

Results

SPME conditions for the analysis of antipsychotic drugs in biofluids were optimized in order to select the most effective extraction phase, and finally C18 blades were chosen, as they provided good recoveries and the lowest carryover effect. Desorption of the analytes was performed in optimal desorption mixture

(ACN/H₂O/FA, 80:20:0.1, v/v/v) and reached equilibrium within only 10 min. The obtained kinetic profile of antipsychotic drugs in saliva samples was compared to the kinetic profiles of these drugs obtained for serum samples. Although the level of antipsychotics in saliva samples was much lower (ng/mL) in comparison to the concentration of drugs measured in serum samples (µg/mL), the kinetic profile of analyzed compounds was similar. Therefore, SPME may be implemented for direct monitoring of the level of antipsychotic drugs in non-invasively collected saliva samples and can be used as alternative to traditional methods utilizing serum samples for the measurements of these drugs.

Conclusions

High-throughput method based on SPME technique coupled to LC-MS/MS was developed for quantitative analysis of five antipsychotics in human saliva and serum samples. The optimized SPME method facilitates non-invasive monitoring of analyzed drugs present at low level in saliva samples, without the need for blood collection. SPME allows for a control of drug intake and dose modification, and may be a part of personalized therapy during long-term treatment.

Novel Aspect

SPME technique for direct monitoring of therapeutic level of atypical antipsychotic drugs in saliva and serum samples was developed.

References

1. Berling I., Buckley N. A., Isbister G., Br J ClinPharmacol. 82, 249-54 (2016).
2. Bessonneau V., Boyaci E., Maciazek-Jurczyk M., Pawliszyn J., Anal ChimActa. 856, 35-45 (2015).
3. Roszkowska A., Tascon M., Bojko B., Goryński K., Dos Santos P. R., Cypel M., Pawliszyn J., Talanta 183, 304-310 (2018).

TP-33 / MULTI-ATTRIBUTE MONITORING (MAM) TO IDENTIFY DIFFERENCES IN TRASTUZUMAB FROM 2 MANUFACTURERS

Ferran Sanchez⁽¹⁾ - ***Chris Nortcliffe***⁽²⁾ - ***Kerstin Pohl***⁽³⁾ - ***Jörg Dojahn***⁽³⁾ - ***Sibylle Heidelberger***⁽²⁾

SCIEX, Biologics, Barcelona (1) - SCIEX, Biologics, Warrington (2) - SCIEX, Biologics, Darmstadt (3)

Keywords: Mass Spectrometry, Biologics, Trastuzumab, MAM, Biosimilars

Introduction

Biosimilars are a growing industry, which require comparison of the safety and efficacy of the new biotherapeutic to the original. Certain CQA's, such as high mannose content glycans, can have a detrimental effect on the function of the product. Using trastuzumab from two different manufacturers, a comparison of the biosimilars was performed, acquiring data using SWATH® and processing using BioPharmaView™ MAM workflow to identify the levels of change in the sample.

Methods

Peptide digests of trastuzumab from two manufacturers, were run using an ExionLC™ AD with a Phenomenex® Kinetex® column C18, 100A, 1.7 µm, 2.1 x 150 mm at 60°C. A: 0.1 % formic acid in water, B: 0.1% formic acid in acetonitrile using a linear gradient of 5 to 40 over 120 min at a flow rate of 0.2 mL/min using SWATH acquisition. Trastuzumab (intact and reduced) were acquired using a 20 min gradient. All data was processed using BioPharmaView MAM and intact workflow.

Results

Trastuzumab peptide mapping was acquired on the X500B using SWATH acquisition. The data was processed in BioPharmaView with the glycosylations setup as attributes for monitoring and the subsequent

glycosylation percentages were calculated and compared against each other for the two manufacturers versions of trastuzumab. Differences in the glycosylations were identified at the peptide level. The same samples were acquired as intact protein and the levels of glycosylations and the variances were compared to the data acquired in peptide mapping. The result was a confirmation that some levels of glycans differed between the two samples showing how the workflow can be used to monitor for changes in biologic in production or can be used to compare biosimilar to the original biotherapeutic.

Conclusions

With the increasing number of biosimilars under development; robust, reliable workflows are required to fully characterize any Critical Quality Attributes which may differ from the originator. Certain process-specific post-translational may have a negative effect on the efficacy and clearance of the therapeutic therefore a Multiple Attribute Method allows monitoring of these in a single workflow.

Novel Aspect

Comparison of 2 biosimilar trastuzumab using multi-attribute monitoring.

TP-34 / LC-MS/MS ANALYSIS AND RAT PHARMACOKINETICS OF THE NOVEL CYCLIN-DEPENDENT KINASE INHIBITORS BP-14 AND BP-20

Jitka Šíroká⁽¹⁾ - *Martina Čečková*⁽²⁾ - *Lubor Urbánek*⁽¹⁾ - *Vladimír Kryštof*⁽¹⁾ - *Tomáš Gucký*⁽¹⁾ - *Jakub Hofman*⁽²⁾ - *Miroslav Strnad*⁽¹⁾ - *František Štaud*⁽²⁾

Laboratory of Growth Regulators, Palacký University, Olomouc (1) - Faculty of Pharmacy in Hradec Králové, Charles University, Hradec Králové (2)

Keywords: cyclin-dependent kinase inhibitor; BP-14; BP-20; pharmacokinetics; tissue distribution

Introduction

BP-14 and BP-20 are novel cyclin-dependent kinase inhibitors, structurally related to roscovitine, with potent antiproliferative activity [1]. A simple LC-MS/MS method for determining them in biological material was developed and validated following EMA guidelines. The method was applied in a pharmacokinetic study in rats following i.v. and i.d. administration including plasma pharmacokinetics, tissue distribution and excretion (renal and biliary).

Methods

Chromatographic separation was performed on reversed phase BEH C18 column (100 x 2.1 mm, 1.7 μm) by gradient elution with mobile phases composed of 15 mM ammonium formate pH 4.0 and methanol at flow rate 0.25 mL/min at 40°C. The analytes were detected based on their multiple reaction monitoring transitions in positive electrospray ionisation mode m/z 473.07 > 157.93 for BP-14, m/z 499.62 > 184.2 for BP-20 and m/z 355.5 > 90.86 for roscovitine (internal standard).

Results

The method provided good linearity within concentrations 1-10000 nmol/L for BP-14 and 10-25000 nmol/L for BP-20 in plasma and 100-10000 nmol/L for BP-14 and 10-25000 nmol/L for BP-20 in bile and urine ($R^2 > 0.9989$). The results of validation fit within the EMA acceptance limits in terms of within-run and between-run accuracy and precision. Both compounds were stable in a set of stability tests in matrices (plasma, bile and urine).

Plasma concentrations of both compounds after i.v. administration followed a typical bi-exponential decline, consisting of a rapid distribution and slow elimination phase. The bioavailability after i.d. administration was determined 0.63% for BP-14 and 1.58% for BP-20. The distribution into tissues (adipose tissue, kidney, lungs, spleen, liver, muscle, brain) was the highest for the adipose tissue and the lowest for

brain. The biliary excretion of the parent BP-14 and BP-20 accounted for 4.81% and 10.6% of the doses, respectively, and renal excretion for less than 0.5% in both cases.

Conclusions

Similar pharmacokinetic profiles of both compounds were obtained after i.v. administration, with a rapid initial decrease in plasma concentrations, intensive distribution to adipose tissue and approximately 5-10% excretion of the parent drugs through urine and bile. Their bioavailability after i.d. administration was low.

Novel Aspect

The results represent first information on pharmacokinetics of BP-14 and BP-20 in rats and will help to design a new generation of drugs with strong anticancer activity.

References

1. Gucky, T., Jorda, R., Zatloukal, M., Bazgier, V., Berka, K., Reznickova, E., Beres, T., Strnad, M., Krystof, V., *Journal of Medicinal Chemistry*, 56, 6234-6247 (2013)

TP-35 / IDENTIFICATION OF 5-FLUOROCYTOSINE AND ITS METABOLITES IN COMPLEX SAMPLES BY LIQUID CHROMATOGRAPHY COUPLED WITH MASS SPECTROMETRY

Pavel Škvára⁽¹⁾ - *Miroslava Matúšková*⁽²⁾ - *Erika Ďuriníková*⁽²⁾ - *Andrea Vojs Staňová*⁽¹⁾

Department of Analytical Chemistry, Faculty of Natural Sciences, Comenius University in Bratislava, Bratislava, Slovak Republic (1) - Institute of Experimental Oncology, BMC, Slovak Academy of Sciences, Bratislava, Slovak Republic (2)

Keywords: mass spectrometry, gene therapy, chemotherapeutics, metabolites, complex samples

Introduction

5-fluorocytosine (5FC) is an anti-mycoticum which can be used as precursor for anti-cancer drug 5-fluorouracil (5FU) in cancer treatment by the procedure "suicide gene therapy" [1]. In recent years, the use of high performance liquid chromatography (HPLC) coupled with mass spectrometry (MS) is dominated for purpose of identification and determination of 5FC and its metabolites in complex samples [2], [3].

Methods

All results were obtained by high resolution LC-MS-IT-TOF mass spectrometer (Shimadzu, Kyoto, Japan) with electrospray ionization. Chromatographic separation was performed on Ascentis C18 (100x2.1 mm; 3 µm) column (Sigma-Aldrich). The samples for analysis was prepared and delivered by the Laboratory of Molecular Oncology.

Results

The developed HPLC-ESI-IT-TOF-MS method was applied to the analysis of samples of mesenchymal stem cells, tumor cells and their control samples to which 5FC was not added. In spectra obtained by HPLC-ESI-IT-TOF-MS analysis of samples (media of tumor or mesenchymal stem cells) were present molecular ions which corresponding to the 5FC in positive ionization mode and to the 5FU and 5-fluoro-2-deoxyuridine monophosphate (5FdUMP) in negative ionization mode. The presence of 5FC's metabolites was conditional upon to presence genes which encoding excretion of enzymes cytosine deaminase (CD) or uracil phosphoribosyltransferase (UPRT). In spectra obtained from analysis of sample without addition of gene CD or CD/UPRT were present only molecular ions corresponding to 5FC and none molecular ions which corresponding to its metabolites in negative ionization mode.

Conclusions

The paper dealt with the possibilities of high performance liquid chromatography coupled with mass

spectrometry for identification of 5-fluorocytosine and its major metabolites in biological samples. Developed HPLC-ESI-IT-TOF-MS method was successfully applied for analysis of mesenchymal stem cells, tumor cells and their control samples.

Novel Aspect

Obtained results can bring new look on the treatment efficiency of metabolic conversion 5FC during suicide gene therapy.

References

1. Zhang J., Kale V., Chen M., An Official Journal of the American Association of Pharmaceutical Scientists, 17, 102–110 (2015)
2. Ishii H., Shimada M., Yamaguchi H., Mano N., Biomedical Chromatography, 30, 1882–1886 (2016)
3. Serve K. M., Yáñez J. A., Remsberg C. M., Davies N. M., Black M. E., Biomedical Chromatography, 24, 556–561 (2010)

Acknowledgements

The research was supported by the grant of the Slovak Research and Development Agency (APVV-16-0124).

TP-36 / EFFICIENT DEREPLICATION WITH HIGH RESOLUTION HIGH ACCURACY MASS SPECTROMETRY FOR THE DISCOVERY OF NOVEL BIOACTIVE NATURAL PRODUCTS

Lijiang Song

Department of Chemistry, University of Warwick, Coventry

Keywords: High resolution, mass spectrometry, dereplication, drug discovery

Introduction

Antimicrobial resistance poses a serious threat worldwide and new antibiotics are urgently required. Meanwhile, the pace of discovery of novel antibiotics has slowed down significantly recently due to the frequent re-discovery of known compounds. Microbial genome mining as a new approach can potentially rejuvenate the stalled antibiotics discovery pipeline. High resolution mass spectrometry provides one of the most efficient ways to carry out dereplication.

Methods

Secondary metabolites produced by producing microorganisms are extracted from spent fermentation medium and analyzed with high resolution LC-MS/MS. Molecular formulae for all major peaks identified in chromatograms are generated based on isotopic pattern and accurate measurement, these will be used for database searching, both known and potentially novel compounds can be further identified through MS/MS and UV absorbance.

Results

Rapid advances in bioinformatics enables the complete or partial structure of novel natural product to be predicted based on gene sequence, combining with genetic manipulation of bacterial strains, the final biosynthetic products can be identified through a LC-MS/MS based comparative metabolomics approach. The discovery and biosynthetic investigation of coelimycin P1, an unusual yellow-pigmented metabolic product of the cpk cryptic polyketide biosynthetic gene cluster from *Streptomyces coelicolor* M145 will be discussed.

Novel Aspect

Efficient de-replication with high resolution high accuracy Q-TOF mass spectrometry is critical for future bioactive natural product discovery.

References

1. Brauner A., Fridman O., Gefen O., Balaban N., Nat. Rev. Microbiol., 2016, 14, 320-330
2. Laureti L., Song L., Huang S., Corre C., Leblond P., Challis G. L. and Aigle B. . Proc. Nat. Acad. Sci. USA. 2011, 108, 6258-6263
3. Gomez-Escribano J, Song L., Fox D., Bibb M. J., Challis G. L., Chem. Sci, 2012, 3, 2716-2720

TP-37 / COLLISION-INDUCED DISSOCIATION STUDY OF ANTIMICROBIAL PEPTIDE DENDRIMERS

Aura Tintaru - Laurence Charles

AIX-Marseille Université - CNRS, Institut de Chimie Radicalaire UMR7273, Marseille

Keywords: Antimicrobial Peptide dendrimers, Collision-induced dissociation, MS/MS sequencing

Introduction

Peptide dendrimers are non-natural branched peptides with multibranching topologies that exhibit interesting anti-microbial properties [1]. In particular, dendritic architectures composed of lysine (K) and leucine (L) were demonstrated to be active against a broad panel of multidrug-resistant pathogens [2]. The renowned efficiency of CID for MS/MS sequencing of linear peptides was evaluated for structural characterization of a KL peptide dendrimer.

Methods

High resolution MS and MS/MS experiments were performed after positive ion mode electrospray ionization (ESI), using a Synapt G2 HDMS instrument. Lysine-Leucine (KL) peptide dendrimers were dissolved in acidified methanol to promote formation of protonated species.

Results

MS/MS was previously demonstrated to be very powerful at characterizing dendritic structures [3], shown to exhibit a typical dissociation behavior consisting of iterative losses of repeating building blocks from outside branches and then internal arms down to the inside core [4]. MS/MS data recorded on the KL-peptide dendrimer suggested a different scenario, with the dissociating precursor ion to be considered as a peptidic backbone with pendant ramified segments rather than a dendritic structure. On the one hand, CID data recorded at low collision energy allowed the main backbone to be sequenced, and hence the branching points to be localized. On the other hand, raising activation energy permitted to characterize the structure of the pending peptidic segments. Moreover, charge states observed for product ions allowed some insights in the actual location of protons in the quadruply charged studied species.

Conclusions

This study showed that KL-peptide dendrimers behave as a branched peptide rather than a dendritic structure upon CID. Dissociation rules to be established from MS/MS data obtained for this first generation molecule are currently evaluated for higher generation KL peptide dendrimers.

Novel Aspect

First report of MS/MS structural characterization KL-peptide dendrimers.

References

1. Mintzer M. A., Dane E. L., O'Toole G. A., Grinstaff M. W., Mol. Pharmaceutics 9, 342 (2012).
2. Stach M., Siriwardena T. N., Köhler T., van Delden C., Darbre T., Reymond J. -L., Angew Chem Int Ed Engl., 53, 12827 (2014).
3. Tintaru A., Chendo C., Monnier V., Bouillon C., Quelever G., Peng L., Charles L., Int. J. Mass Spectrom. 308, 56 (2011).

4. Tintaru A., Ungaro R., Liu X., Chen C., Giordano L., Peng L., Charles L., Anal. Chim. Acta 853, 451 (2015).

TP-38 / MASS SPECTROMETRY CHARACTERIZATION OF AN NMS PROPRIETARY DRUG CONJUGATED ANTIBODY (ADC)

Sonia Troiani⁽¹⁾ - **Rita Perego**⁽¹⁾ - **Ulisse Cucchi**⁽¹⁾ - **Paolo Orsini**⁽²⁾ - **Matteo Salsa**⁽²⁾ - **Michele Caruso**⁽²⁾ - **Barbara Valsasina**⁽¹⁾

Nerviano Medical Sciences srl, Oncology/Biotechnology, Nerviano (MI) (1) - Nerviano Medical Sciences srl, Oncology/Medicinal Chemistry, Nerviano (MI) (2)

Keywords: Antibody-Drug conjugates (ADCs), Cysteine-conjugation, Hydrophobic Interaction Chromatography (HIC), Mass Spectrometry

Introduction

Antibody-Drug Conjugates (ADCs) combine the target specificity of monoclonal antibodies to the high cytotoxicity of a small molecule toxin selectively killing tumor cells while minimizing toxicity to normal tissues. ADCs are generated by covalently binding the selected toxin through a linker to specific antibody residues, usually Cysteine or Lysine, which increases the complexity of the products generating heterogeneous mixtures of isomers and isoforms

Methods

The characterization of the ADC product is fundamental because drug loading and distribution can affect the safety and efficacy of the ADC. Mass spectrometry can be applied for the identification of the ADC modified residues and, in combination with chromatographic techniques as UV-LC or HIC, for the determination of the Drug-Antibody Ratio (DAR) and the distribution of the toxin in the different populations generated by the cysteine conjugation process

Results

At NMS a new proprietary molecule suitable for conjugation has been developed and used to derivatize a commercial antibody to generate a cysteine conjugated ADC with an average DAR of ~3.5. The addition of the drug at the interchain cysteine residues results in a heterogeneous population of ADCs that differ in the number of drugs per antibody and in the presence of positional isomers.

The peptide mass fingerprinting strategy has been applied for the identification of the modified residues after the complete enzymatic digestion of the ADC and to characterize the ADC at the peptide level looking at primary sequence and post-translational modifications.

The DAR distribution analysis and the identification of the positional isomers present in the sample have been achieved by the combination of HIC and RPLC-MS. The intact DAR species separated by HIC on the bases of their hydrophobicity have been collected and analysed in denaturing conditions by RPLC-MS in order to separate and identify the ADC fragments composing each isomeric form

Conclusions

The application of analytical technologies, especially mass spectrometry, have been applied for the analysis of a cysteine conjugated ADC produced at NMS. These technologies represent very useful tools for the characterization of these heterogeneous products during the drug development process to guarantee the batch to batch reproducibility

Novel Aspect

Mass spectrometry, in combination with allied separation techniques, has been applied for the characterization of an NMS proprietary drug conjugated antibody

TP-39 / ASSESSING THE DRUGABILITY OF PTP1B NON-CATALYTIC SITES WITH MASS SPECTROMETRY

Ana Villar Garea⁽¹⁾ - *Christian Buning*⁽¹⁾ - *Joerg Czech*⁽¹⁾ - *Matthias Dreyer*⁽¹⁾ - *Gerhard Hessler*⁽¹⁾ - *Stefan Matthias Petry*⁽¹⁾ - *Norbert Tennagels*⁽²⁾ - *Stefan Welte*⁽²⁾

Sanofi-Aventis Deutschland GmbH, Integrated Drug Discovery, Frankfurt am Main (1) - Sanofi-Aventis Deutschland GmbH, Diabetes Research, Frankfurt am Main (2)

Keywords: inhibitor, drug discovery, target engagement, cysteine alkylation, phosphatase

Introduction

The identification of therapeutic inhibitors for protein tyrosine phosphatase 1B (PTP1B) is an attractive approach to treat type II diabetes, but it has failed so far due to insufficient selectivity [1]. Covalent inhibitors that target residues only present in PTP1B might overcome these limitations. We studied the reactivity of cysteines in PTP1B and designed and characterized potential ligands.

Methods

Mass spectrometry: 2-(4'-(iodoacetamido)anilino)naphthalene-6-sulfonic acid (IAANS) was used for modification of accessible cysteine residues of PTP1B. The number and position of the reactive residues were determined by LC-MS1 and LC-MS2. Reaction of PTP1B with the designed ligands was studied in the same way. Other methods: molecular modelling, in vitro phosphatase assay, X-ray crystallography.

Results

First, the accessibility of different cysteine residues was determined by incubating PTP1B with IAANS under native conditions and analyzing the resulting products by LC-MS. We found Cys32 as potential target for specific inhibition of PTP1B. Next, a library of molecules was designed for covalent inhibition of PTP1B and was screened with an in vitro enzymatic assay. Several inhibitors were selected for mass spectrometry experiments to evaluate their reactivity towards Cys32. Despite a having a high potency, no alkylation of Cys32 with these ligands could be detected under the same conditions as used for IAANS. X-ray crystallography analyses revealed that, although the inhibitors occupy the intended pocket, their orientation hinders the reaction with Cys32.

Conclusions

We identified Cys 32 as potential target for PTP1B inhibition and designed inhibitors to alkylate it. Although these substances are very potent in enzymatic in vitro assays and bind the expected pocket, neither X-ray crystallography nor mass spectrometry could demonstrate covalent modification of PTP1B. Therefore we conclude that the observed activity is due to non-covalent interaction.

Novel Aspect

Use of IAANS to test the accessibility of cysteine residues; Handling a native protein (6 reduced Cys) without reducing agent

References

1. Popov D. Biochemical and Biophysical Research Communications 410, 377–381 (2011)

TP-40 / ULTRAFAST DETECTION OF DRUGS AND METABOLITES IN URINE BY FLOW INJECTION ANALYSIS (FIA) COUPLED TO MAGNETIC RESONANCE MASS SPECTROMETRY (MRMS)

Matthias Witt⁽¹⁾ - *Markus Godejohann*⁽²⁾ - *Aiko Barsch*⁽³⁾

Bruker Daltonik GmbH, MRMS Solutions, Bremen (1) - Bruker BioSpin GmbH, Hyphenation AIC, Rheinstetten (2) - Bruker Daltonik GmbH, Metabolomics, Bremen (3)

Keywords: Metabolomics, Drug detection, Flow injection analysis, Screening

Introduction

Detection of metabolites and drugs in body fluids such as plasma or urine by LCMS is a routine method in metabolomics and doping analysis. Routine UPLC-MS measurements are performed typically in 15 min. Therefore, the number of analyzed samples is highly limited. In this work, a fast method for detection of drugs and their metabolites in urine using FIA-MRMS is presented. Roughly 250 samples can be measured in 24h using this technique.

Methods

Six pooled urine samples were purified by SPE using Merck LiChrolutEN SPE cartridges. Samples were extracted with methanol from SPE cartridges and diluted 1 to 100 for FIA. Each samples was analyzed in 5 minutes by FIA-MRMS using a solarix 2xR (Bruker Daltonik, Bremen) in ESI using positive and negative ion mode. Analysis of data was performed with MetaboScape 3.0 (Bruker Daltonik, Bremen).

Results

The data of the ESI(+) and ESI(-) were combined for feature analysis. More than 2100 features have been found for the pooled urine samples. More than 90% of the detected features could be assigned with a molecular formula. 300 drug candidates have been detected in the urine samples using a HMDB urine database with a mass error tolerance of only 0.5 ppm. The detected drugs have been compared with the medication of the patients. Several drugs have been found only in one or a few pooled urine samples. By comparing the relative abundances of features of all samples, possible metabolites of drugs could be identified.

Conclusions

Drugs and their metabolites can be detected by FIA-MRMS in a few minutes. This workflow is much fast than the conventional workflow using UPLC-MS. This method could even be used for quantification when internal drug standards are added. Due to the complexity of the samples ultra-high mass resolution as well as very accurate mass detection is a prerequisite for this workflow.

Novel Aspect

FIA-MRMS can be used for fast detection of drugs and metabolites in urine.

TP-41 / EFFICIENT PHOSPHOLIPIDS REMOVING CLEAN-UP SPE-LC-MS/MS METHOD FOR SIMULTANEOUS PLASMA QUANTIFICATION OF 11 TYROSINE KINASE INHIBITORS AND CAFFEINE. APPLICATION TO THERAPEUTIC DRUG MONITORING

Dóra Koller^(1,2) - *Viktoryia Vaitsekhovich*⁽²⁾ - *Cecile Mba*⁽²⁾ - *Juan Luis Steegmann Olmedillas*⁽³⁾ - *Francisco Abad-Santos*^(1,2) - *Aneta Wojnicz*^(1,2)

Clinical Pharmacology Department, Hospital Universitario de la Princesa, Instituto Teófilo Hernando, Universidad Autónoma de Madrid, Instituto de Investigación Sanitaria la Princesa (IP) (1) - Fundación de la Investigación Biomédica, Instituto de Investigación Sanitaria la Princesa (IP) (2) - Hematology Department, Universidad Autónoma de Madrid, Instituto de Investigación Sanitaria la Princesa (IP), Madrid, Spain (3)

Keywords: liquid chromatography-tandem mass spectrometry, tyrosine kinase inhibitors, μ -elution solid phase extraction, endogenous phospholipids, therapeutic drug monitoring

Introduction

Therapeutic drug monitoring (TDM) can improve efficacy and safety (1). Thus, we developed a liquid chromatography-tandem mass spectrometry (LC-MS/MS) method for monitoring 11 tyrosine kinase

inhibitors (TKIs) used in chronic myeloid leukemia (imatinib, dasatinib, nilotinib, bosutinib, ponatinib), polycythemia vera (ruxolitinib), chronic lymphocytic leukemia (ibrutinib) and rheumatoid arthritis (filgotinib, tofacitinib, baricitinib, peficitinib).

Methods

Only 200 µL of human plasma was needed for the extraction of analytes and their stable isotope-labeled internal standards. µ-elution solid phase extraction was optimized and compared to simple protein precipitation. Gradient elution on Poroshell 120 EC-C18 column, flow rate of 0.6 mL/min and 60°C was applied for compounds separation. The mobile phase consisted of 0.1% formic acid in water, pH=2 (solution A) and 0.1% formic acid in acetonitrile (solution B).

Results

Acquisition lasted for 8 minutes followed by a re-equilibration time of 4 minutes. Mass spectrometer was operating in the positive ionization mode. Dynamic multiple reactions monitoring scan was applied to improve method sensitivity. Endogenous phospholipids can strongly affect MS analysis. Therefore, the monitoring of m/z 184 > 184 (glycerophosphocholines) and m/z 104 > 104 (lysophosphatidylcholines) as common in-source collision-induced dissociation ion fragments, were included in the method. Caffeine was added to the assay as well. According to the recommendations of regulatory agencies [2,3,4], tests of precision, accuracy, recovery, matrix effect, process efficiency, stability, sensitivity (great LLOQs), and selectivity were performed and enabled us to validate the method. Sample preparation applied in the present assay showed high efficiency in phospholipids removal. More than 99% of main plasma lysophosphatidylcholines and more than 89% of glycerophosphocholines were eliminated compared to protein precipitation.

Conclusions

This method enables simultaneous plasma monitoring of 11 TKIs, caffeine and two main plasma phospholipids. High effectiveness in endogenous phospholipids elimination ensures the reliability of the assay. Present LC-MS/MS method is currently used in our clinical practice, being applied to TDM of dasatinib, imatinib, nilotinib and ponatinib. TDM of TKIs helps to individualize dose adjustment and manage adverse effects in patients.

Novel Aspect

The method monitors endogenous phospholipids (glycerophosphocholines and lysophosphatidylcholine) as well as caffeine in addition to 11 TKIs in a simple run.

References

1. Yu H, Steeghs N, Nijenhuis CM, Schellens JH, Beijnen JH, Huitema AD. Clin Pharmacokinet. 53, 305-25 (2014)
2. FDA, Homepage, Guidance for Industry: Bioanalytical Method Validation, 2001, <http://www.fda.gov/cder/guidance/4252fnl.htm>.
3. EMA, Homepage, Guideline on Bioanalytical Method Validation. EMEA/CHMP/EWP/192217/2009,2011, <http://www.ema.europa.eu/docs/enGB/documentlibrary/Scientificguideline/2011/08/WC500109686.pdf>.
4. ICH ICH Harmonized Tripartite Guideline. Validation of Analytical Procedures: Text and Methodology Q2 (R1), 1994, https://www.ich.org/fileadmin/Public_Web_Site/ICH_Products/Guidelines/Quality/Q2_R1/Step4/Q2R1_Guideline.pdf.

TP-42 / QUANTIFICATION OF MAJOR PEANUT PROTEIN ALLERGEN ISOFORM GROUPS FROM PEANUT PROTEIN EXTRACTS

Kim Alving - Pierre DeMontigny - Bing Wang

Sanofi, Analytical Research and Development, Waltham, United States

Keywords: Absolute peanut protein quantification LC/MS/MS

Introduction

Peanut allergy is a potentially fatal condition that has seen a rise of 50% in the recent 15 years currently affecting 0.6-1.3% of the US population. Desensitization by exposure to peanut allergen doses is being evaluated in several clinical trials. Accurate quantification of the peanut allergens is key to successful desensitization. We here report on the development and qualification of a LC/MS based method for absolute quantification of peanut protein content.

Methods

Discovery analysis of peanut protein extract PPE digests was carried out as nanoLC Orbitrap MS/MS analyses, and suitable peptides were selected for quantification. PPE digest allergen quantification was carried out on a triple quadrupole MS in MRM mode. Unlabeled and heavy isotope labeled peptides (internal standards) were used to generate calibration curves. Intact peanut protein standards were applied for accurate determination of total process recovery.

Results

Highlights of the peanut proteome discovery efforts and quantification method development rationale, selection of representative peptides, development and qualification results for the analytical methods to quantify four major allergen proteins (Ara h1, h2, h3 and h6) including linearity, accuracy and reproducibility (triple replicate, triple inter day replicate, as well as reproducibility), will be presented together with the results from the analysis of several PPE research batches. The average recovery was determined to be 61% applying intact allergen standard protein. The intra- and inter-day CV% was below 6.6%, and accuracy is within +/- 10% . The stability data for freeze-thaw and 2-hour room temperature stability and 12 hour protein extract digest stability were all below 15%.

Conclusions

An LC/MS/MS based method for quantification of four major peanut allergens from peanut extracts has been developed and qualified. The method was found to have good sensitivity, reproducibility, and dynamic range and thus to be suitable for monitoring the concentrations of the major peanut allergen proteins in research batches obtained from peanut extracts. The method was designed to detect and quantify the isoforms of each of the 4 major peanut protein allergens.

Novel Aspect

Absolute quantification of PPE with multiple isoforms by LC/MS/MS. Application of intact proteins for accurate recovery determination

TP-43 / VALIDATION OF BREATH BIOPSY PLATFORM FOR LONGITUDINAL MEASUREMENT OF VOCS IN BREATH USING TD-GC-MS

Jasper Boschmans - Rob Smith - Marc van der Schee - Simon Kitchen - Billy Boyle - Max Allsworth

Owlstone Medical Ltd., Cambridge, United Kingdom

Keywords: Breath Research, TD-GC-MS, Volatile Organic Compounds, Pharmacokinetics, Breath Biopsy

Introduction

Exhaled breath contains thousands of Volatile Organic Compounds (VOCs) which are products of metabolic activity, and promising biomarkers for a range of diseases. Breath also includes VOCs produced during the metabolism of pharmaceuticals and other xenobiotics.

Evaluating longitudinal changes in the breath VOC profile is relevant for a wide range of applications such as disease monitoring, measuring response to therapeutic interventions, assessing effects of environmental exposures, and studying pharmacokinetics. This study describes how the Breath Biopsy platform can be used to capture multiple breath samples over time, uncovering detailed changes in the concentration of volatile metabolites present in breath.

We present two main experiments: first we analyze VOCs in breath following ingestion of a peppermint capsule and show that Breath Biopsy can be used to observe the “washout curve” in for target compounds in a single individual using repeated, robust breath collection and analysis over a period of 8 hours. Secondly, we examine the biological variation in a longitudinal study where the washout experiment is repeated in the same individual multiple times over the course of 5 weeks.

Methods

Breath samples were collected using the ReCIVA Breath Sampler and analyzed using the Breath Biopsy platform in the Breath Biopsy Clinical Laboratory (Owlstone Medical Ltd, UK). Samples were pre-purged to remove excess water and desorbed using a TD100-xr thermal desorption autosampler (Markes International) and transferred onto a VF-5ms column (60m x 0.25mm x 0.25 μ m; Agilent Technologies) using 1:2 split injection. Chromatographic separation was achieved via a programmed method (40-310°C in 60 min. at 1.8 mL/min.) on a Trace 1310 GC oven (Thermo Fisher Scientific) and mass spectral data acquired using an electron ionization time-of-flight (*i.e.* EI-TOF) BenchTOF HD mass spectrometer (also Markes International). Raw Markes data files were converted from using TOF-DS (Markes International) and MassHunter Quant (Agilent Technologies) was used for further data processing.

Results

Analysis of breath captured 30 minutes after consumption of the peppermint capsule shows a large increase in the VOCs α -pinene, β -pinene, limonene, eucalyptol and p-menthan-3-one compared to baseline pre-ingestion controls. Breath collections made every 30 minutes after this initial capture show a washout curve for all peppermint-related compounds, decreasing to baseline levels. A wide range of fold changes was observed across all peppermint-related compounds, with most abundant GC peaks not necessarily resulting in highest fold changes. In contrast, breath metabolites *e.g.* acetone that are not peppermint-related show little or no change compared to control.

From the washout experiment, three time points were selected (pre-ingestion control, peak at 45 minutes and plateau at 3 hours) for inclusion in a longitudinal study where the washout experiment was repeated multiple times over 5 weeks.

In the longitudinal study, fold changes for all compounds are presented for the peak and plateau time points, relative to the corresponding pre-ingestion control. Breath metabolites acetone and isoprene show only small differences between peak and plateau fold change as expected, however large fold changes are observed for the peppermint-related VOCs such as α -pinene and β -pinene, although some variation is observed for peak data points between washout experiments.

Conclusions

The study demonstrates that the Breath Biopsy platform can be used to study longitudinal changes of exhaled VOCs in a reliable and reproducible way. The peppermint washout experiment is a useful way to assess platform performance and future work will investigate the biological variability that exists for uptake and breakdown of substrate. The Breath Biopsy platform unlocks potential new use cases for breath analysis in metabolomics and pharmacokinetics studies.

Novel Aspect

Breath-based platform allowing standardized sample collection and pre-concentration using the ReCIVA Breath Sampler, followed by analyte separation and data acquisition via a TD-GC-MS methodology.

References

1. M. M. L. Steeghs, *Int. J. Mass. Spectrom.*, 263, 204-212 (2007).

2. L. Schmidt, Arch. Toxicol., 91, 677-687 (2017).
3. L. M. Heaney, Doctoral Thesis submitted in partial fulfilment of the requirements for the award of Doctor of Philosophy of Loughborough University (UK), (2016).
4. L. M. Heaney, Bioanalysis, 8, 1325-1336 (2016).
5. J. Beauchamp, J. Breath Res., 4, 026006 (2010).

TP-44 / ELIMINATION PROFILE OF BETAMETHASONE AND METABOLITES IN URINE AFTER DIFFERENT ADMINISTRATION ROUTES

Sergi Coll⁽¹⁾ - **Xavier Matabosch**⁽¹⁾ - **Nuria Monfort**⁽¹⁾ - **Elida Alechaga**⁽¹⁾ - **Clara Pérez-Mañá**⁽²⁾ - **Julián A Mateus**⁽²⁾ - **Jordi Monfort**⁽³⁾ - **Jone Llorente-Onaindia**⁽³⁾ - **Rosa Ventura**⁽¹⁾

Institut del Mar d'Investigacions Mèdiques, Catalanian anti-doping laboratory, Barcelona, Spain (1) - Institut del Mar d'Investigacions Mèdiques, Human Pharmacology and Clinical Neurosciences Research Group, Barcelona, Spain (2) - Institut del Mar d'Investigacions Mèdiques, Cell Research on Inflammation and Cartilage, Rheumatology department, Barcelona, Spain (3)

Keywords: glucocorticoids, anti-doping, administration routes, LC-MS/MS

Introduction

Betamethasone (BET) is a glucocorticoid widely used in sports medicine for the treatment of asthma and acute injuries. BET is prohibited in sport competitions when administered by systemic routes (*e.g.*, oral, intramuscular), and is allowed using local treatments. The aim of the present work was to examine the urinary profile of BET and its metabolites after different administration routes.

Methods

A validated method was used to determine urinary levels of BET, 6 β -hydroxy-BET, and fifteen additional metabolites. Urine samples were hydrolyzed with β -glucuronidase, extracted with ethyl acetate and analyzed by LC-MS/MS. The method was applied to urines collected in clinical studies where BET was administered by topical, intranasal, intra-articular, peri-articular, oral and intramuscular routes to male and female volunteers (n=40).

Results

After topical (10 mg/day for 5 days) and intranasal (320 μ g/day for 3 days) uses, maximal BET concentrations were 6 and 32 ng/mL, respectively. In most of the samples collected after intranasal treatment, BET concentrations were below 30 ng/mL.

Intra-articular and peri-articular single administrations led to maximal BET concentrations of 1420 and 981 ng/mL, respectively. All samples collected 72h after both administrations gave BET concentrations lower than 30 ng/mL.

After a single oral dose of 0.5 mg, the maximal BET concentration was 82 ng/mL, and all samples collected 24h after intake resulted in BET concentrations below 30 ng/mL. Single intramuscular doses of 6 and 12 mg gave maximal BET concentrations of 1446 ng/mL. In most of the samples collected 72h after intramuscular administration, BET concentrations were lower than 30 ng/mL. After both administrations, BET excretion rates reached maximal values during the first hours after intake and gradually decreased.

The excretion profiles of most of the metabolites were similar to those of BET.

Conclusions

As expected, different excretion profiles were obtained depending on the administration route. It will be difficult to distinguish intra-articular and peri-articular uses from intramuscular treatments because excretion profiles obtained after these three administrations routes are similar. Topical and intranasal treatments can be differentiated from oral and intramuscular administrations, using a reporting level for BET.

Novel Aspect

Since BET can be administered by different routes, this information could be useful for the anti-doping community to know how it behaves depending on the administration route.

TP-45 / VARIABILITY AND VALUES IN GLYCOSYLATION ANALYSES OF MONOCLONAL ANTIBODIES: AN INTERLABORATORY STUDY

M. Lorna A. De Leoz - David L. Duewer - Stephen E. Stein

National Institute of Standards and Technology (NIST), Material Measurement Laboratory, Gaithersburg, United States

Keywords: glycomics, glycoproteomics, NISTmAb, glycan, glycopeptide

Introduction

Glycosylation is a critical post-translational modification in monoclonal antibodies. However, it is inherently heterogeneous and challenging to analyze, leading to a proliferation of analysis methods. An interlaboratory study was coordinated by NIST to determine variability in measuring *N*-glycans across laboratories and to derive consensus values for NISTmAb, an IgGk1 reference material. This work describes results from 103 reports worldwide [1].

Methods

Participants performed glycosylation analysis of NISTmAb and mod-NISTmAb using their own method. Labs recorded, in a provided template, methods used and percent fractional abundances of glycans in both samples. A pilot study was conducted to polish the study's protocol. Consensus values and ratios were derived for glycans reported by seven or more labs. Data was analyzed using robust techniques to assess measurement reproducibility and characterize distributions.

Results

A total of 103 reports were submitted by 76 participants from industry, university, research, and government laboratories in Europe, North America, Asia, and Australia. Various methods were used for glycosylation analysis, including mass spectrometry, fluorescence detection, and capillary electrophoresis. Fifty-seven glycan compositions were reported at least six times; some labs differentiated between isomers. Consensus values were calculated for each of the 57 glycans. The number of glycan compositions identified by labs ranged from 4 to 48, with a median of 22. Most participants reported about the same number of glycans for NISTmAb and mod-NISTmAb. To ascertain the ability of labs to demonstrate comparability of two samples, the mod-NISTmAb/NISTmAb ratios were calculated. The ratios for the four most commonly reported glycans agreed within 8% to 30%. Notably, the extent of agreement between a participant's reported values and the study's consensus values is a strong function of the participant's measurement repeatability.

Conclusions

Determining variability of glycosylation profiles and assigning best values from a larger number of laboratories and range of measurement methods provides a baseline for comparison in the rapidly developing field of glycosylation analysis. We thank all the scientists who participated in this study.

Novel Aspect

Variability and values for the glycosylation of a monoclonal antibody reference material from 100+ reports worldwide

References

1. De Leoz, MLA, Duewer, DL, Stein, SE, NIST Internal Report, 8186 (2017).

TP-46 / APPLICATION OF ACOUSTIC MIST IONIZATION MASS SPECTROMETRY (AMI-MS) FOR THE DEVELOPMENT OF KINASE ASSAYS

Deepika Srivastava⁽¹⁾ - *Phillip DeLand*⁽²⁾ - *Jonathan Wingfield*⁽³⁾ - *Jonathan Tart*⁽¹⁾ - *Paul Clarkson*⁽¹⁾ - *Michael Tonge*⁽³⁾

AstraZeneca, Discovery Biology, Cambridge, United Kingdom (1) - Labcyte Inc., Mass Spectrometry Product Management, San Jose, United States (2) - AstraZeneca, Mechanistic Biology & Profiling, Cambridge, United Kingdom (3)

Keywords: acoustic loading; mass spectrometry; high-throughput screening; kinase; ADP Glo

Introduction

We have implemented a platform utilizing acoustics to load samples into a time-of-flight mass spectrometer capable of analysing more than 100,000 samples per day. In this poster we demonstrate the use of our platform to screen a kinase and compare the data with a classical high-throughput screening technology. Since acoustic mist ionization mass spectrometry (AMI-MS) is a direct infusion system, kinases are a challenge due to the enzyme's requirement for MgCl₂.

Methods

An acoustic transducer generates a nebulised spray from the wells of a 384-well plate located on a moving XY-stage. High voltage is applied to a charging cone above the well, inducing charge separation in the sample. The nebulised spray is drawn through a heated transfer interface into the MS and an ion beam is generated. To utilize the platform with a kinase enzyme, ammonium phosphate at high pH is added to precipitate out the Mg²⁺ ions.

Results

AMI-MS was used to build a kinase enzyme assay to identify inhibitors from a subset of the AZ compound collection. We established a protocol using ammonium phosphate as a stop solution at the end of the assay to remove magnesium chloride and reduce suppression in the mass spectrometer. Having established the kinetic parameters for this target we screened a 60,000 compound subset using both ADPGlo™ kits and AMI-MS endpoints. There was a strong correlation between both assay formats, over 2000 compounds were identified as active in both formats. Surprisingly, there were only small clusters of compounds which appeared active in only one of the two assay technologies (38 unique to AMI-MS and 101 unique to ADP Glo™).

Conclusions

Our data demonstrates that AMI-MS is able to screen kinase targets where there is a requirement for magnesium chloride despite this being a direct infusion MS technology. AMI-MS compared very favourably with traditional kinase screening technologies but proved significantly cheaper as there was no requirement to purchase detection kits.

Novel Aspect

This is the first example of AMI-MS being used to screen a kinase target. The data presented clearly demonstrates that AMI-MS is comparable with existing kinase technologies.

TP-47 / A MASS SPECTROMETRY TOOLBOX FOR THE OPTIMIZATION OF K-RAS G12C COVALENT INHIBITORS

Sophie Raibaud - *Armelle Buzy* - *S. Hamon* - *P. Paul* - *Jean-Claude Guillemot* - *V. Lalleman* - *P. Didier* - *Rosalia Arrebola* - *Francis Duffieux* - *Annick Parent* - *Alexey Rak* - *Gary McCort* - *Laurent Debussche* - *Florence Fassy*

Sanofi, Research Platform & Oncology, Chilly-Mazarin & Vitry sur Seine, France

Keywords: Covalent inhibitor, target engagement, K-RAS

Introduction

Covalent inhibitor of K-RAS G12C to fight lung cancer.

The small GTPase K-RAS is the most prevalent oncogene in cancers. The Gly to Cys mutation at position 12 is frequent in lung cancer. The optimization of a drug, covalently bound to Cys12, would target the oncogene without affecting the wild type protein [1]. Mass spectrometry is well suited to the study of covalent complexes and we have implemented a series of tools to support a K-RAS G12C research program.

Methods

Mass spectrometry tools to study covalent complexes.

Several LCMS/MS-based assays were developed either with recombinant protein or with cellular or tissue lysates: intact mass to rapidly screen or fully characterize the binding kinetics (kinact/Ki) of covalent complexes, peptide mapping to determine the binding site, and targeted quantitative proteomics to measure target engagement in culture cells and tissues.

Results

Optimization of K-RAS G12C covalent binders.

To discover novel covalent inhibitors of K-RAS G12C, several thousands of covalent compounds were screened in an assay measuring the extent of covalent binding on a recombinant protein by intact mass LCMS. Hits were validated by peptide mapping experiments to ascertain that the Cysteine 12 residue was the site of interaction between the compounds and K-RAS G12C. During a chemical optimization program, the newly synthesized compounds were evaluated in the intact mass LCMS assay, as a function of inhibitor concentration at a fixed incubation time (EC₅₀) or at several times (kinact/Ki) [2]. Optimized compounds inhibited the growth of a K-RAS G12C dependent cell line. A targeted LCMS-based assay was implemented to directly and quantitatively determine the engagement of KRAS-G12C in cell lysates [3].

Conclusions

Our LCMS toolbox has proven to be very robust and powerful to support the optimization of novel covalent inhibitors of K-RAS G12C. Although other biochemical or cellular assays were implemented for the research program, the LCMS assays provided the most direct and non-ambiguous results.

Novel Aspect

An LCMS intact mass assay as a primary tool to support chemical optimization in a pharmaceutical discovery program.

References

1. Ostrem JM, Peters U, Sos ML et al. K-Ras(G12C) Nature, 503, 548 (2013)
2. Hansen R, Peters U, Babbar A et al. Nature Struct. And Molec. Biol., doi: 10.1038/s41594-018-0061-5 (2018)
3. Patricelli MP, Janes MR, Li L-S et al. Cancer Disc., 6, 316 (2016)

TP-48 / LC-MS/MS METHOD FOR ANALYSIS OF THERAPEUTIC PEPTIDE AXON PEPTIDE 108 IN AADVAC1 DRUG PRODUCT

Jaroslav Galba⁽¹⁾ - **Juraj Piestansky**⁽²⁾ - **Peter Barath**⁽³⁾ - **Branislav Kovacech**⁽²⁾ - **Andrej Kovac**⁽¹⁾

Axon Neuroscience SE R&D, Proteomics & Metabolomics, Bratislava, Slovakia (Slovak Republic) (1) - Axon Neuroscience SE R&D, Structural Proteomics, Bratislava, Slovakia (Slovak Republic) (2) - Slovak Academy of

Keywords: liquid chromatography, tandem mass spectrometry, Alzheimer's disease, AADVac1, chemical digestion

Introduction

Alzheimer's disease (AD) is a neurodegenerative disorder responsible for more than 70% of all dementia cases. This disease is characterized neurologically by progressive loss of memory, behavior, language, and visuo-spatial skills impairments. AD has become the focus of one of the most intensive investigation in the medicinal field because of worldwide universal phenomenon – population aging and considerable economic impact. Novel therapeutic strategies based on biologicals are currently under investigation recently. One of the innovative AD therapeutic in clinical trials is AADvac1 – immunotherapy designed to elicit an immune response against pathologically modified forms of tau protein, developed by Axon Neuroscience [1,2]. Herein, we present combination of simple acidic chemical cleavage (sample preparation step) with effective and very fast LC-MS/MS method for quantification of immunogen Axon Peptide 108 in AADvac1.

Methods

Chemical cleavage:

A strategy based on chemical cleavage at aspartyl residues using dilute formic acid developed by Li et al. [3] was used in this study. Appropriate cleavage of the therapeutic peptide was obtained when the sample was incubated in 2% formic acid at 108°C for 2 h.

LC-MS/MS method:

The quantitative LC-MS/MS experiments were performed with a Waters ACQUITY UPLC system coupled to the Quattro Premier XE triple quadrupole mass spectrometer. ACQUITY UPLC BEH Shield RP18 column (2.1 mm x 100 mm, 1.7 μm particle size) was used for analysis. Column temperature was set for 30°C. Mobile phase A consisted of 20 mM ammonium formate/0.1% FA in MPW, and mobile phase B was acetonitrile. The flow rate was 0.5 mL min^{-1} and the injection volume was 5 μL . Mass spectra were acquired using positive electrospray ionization mode (ESI+) and selected reaction monitoring (SRM) mode. The following MS conditions were applied. The capillary voltage was 3 kV and the source temperature and desolvation temperature were 120 and 450°C, respectively. The cone gas and desolvation gas flowed at 50 and 600 L/h, respectively. Argon was used as collision gas.

Results

The effective non-enzymatic release step of Peptide 108 based on acid hydrolysis with the use of 2% formic acid was successfully tested and implemented. The main benefit of the chemical cleavage procedure performed with formic acid was high level of compatibility of the cleaved sample with the ESI-MS ionization. Moreover, formic acid was a part of used mobile phase. The developed and optimized LC-MS/MS method was then validated according to the ICH guideline Q2(R1) [4]. Calibration curve was linear within the range of 10 – 300 $\mu\text{g} \cdot \text{mL}^{-1}$. Correlation coefficient was higher than 0.99. The intra- and inter-day precision were 0.92-5.03% and 3.05-10.36%, respectively. The recovery ranged in the interval of 99.73-101.10%. The developed and validated method was successfully applied to determine the amount of Axon Peptide 108 in several clinical batches of AADVac1. Three different batches of AADVac 1 drug product were analyzed and the determined amount of the therapeutic peptide ranged from 74.8 to 106.8 $\mu\text{g}/\text{mL}$.

Conclusions

We developed and tested an effective approach of chemical cleavage of peptide conjugated to the high molecular protein carrier which was connected with very effective, fast and simple LC-MS/MS analysis. Developed and validated LC-MS/MS method was successfully applied and the exact amount of Axon Peptide 108 presented in AADVac 1 was determined in final drug product. The obtained results were in good accordance with the declared amount. This relatively rapid, simple and robust method has a potential to be applied in to the field of biopharmaceuticals analysis.

Novel Aspect

A simple and rapid approach based on combination of non-enzymatic cleavage and LC-MS/MS analysis was presented as a powerful tool for quantitation of therapeutic peptide against Alzheimer's disease in innovative biologic AADvac1. The use of chemical hydrolysis was associated with solubilization of the AADvac 1 adjuvant (alhydrogel) which was beneficial in term of simplicity. Moreover, it was clearly demonstrated that chemical hydrolysis is much more faster and cheaper in comparison to convenient approach based on enzymatic cleavage (*e.g.* asparaginase N).

References

1. Novak M., Kontsekova E., Kovacech B., Zilka N., Protein-based therapy and diagnosis of tau-mediated pathology in Alzheimer's disease, WO2013041962 A1, 2013.
2. Kontsekova E., Zilka N., Kovacech B., Novak P., Novak M., Alzheimer's Research & Therapy, 6, 4(2014).
3. Li A., Sowder R. C., Henderson L. E., Moore S. P., Garfinkel D. J., Fisher R. J., Analytical Chemistry, 73, 5395 (2001).
4. ICH Harmonised Tripartite Guideline Q2(R1), Validation of Analytical Procedures: Text and Methodology, November 2005.

TP-49 / ANALYSIS OF HOST CELL PROTEINS IN ANTIBODY PREPARATIONS USING PASEF

Michael Greig⁽¹⁾ - **Stuart Pengelley**⁽²⁾ - **Guillaume Tremintin**⁽¹⁾ - **Detlev Suckau**⁽²⁾ - **Waltraud Evers**⁽²⁾

Bruker, Daltonics, San Jose, United States (1) - Bruker, Daltonics, Bremen, Germany (2)

Keywords: TimsTof, PASEF, HCP, proteomics, host cell protein

Introduction

In the arena of biotherapeutic analysis, the analysis of host cell proteins (HCP) at the ppm level is critical. ELISA is currently the gold standard for QC applications, but the advantages of mass spectrometry are abundant. We show how PASEF (parallel accumulation and serial fragmentation), as implemented on the timsTOF PRO QTOF, can be applied to HCP analysis for highly sensitive detection with enhanced data quality.

Methods

The NISTmAb Reference Material 8671 and the Universal Proteomics Standard (UPS1, Sigma) were reduced, and digested. Peptides were separated on an Intensity Solo 2 1.8µm C18 column using an Elute UHPLC coupled to a timsTOF Pro ion mobility QTOF mass spectrometer (all Bruker Daltonics). A 150 minute gradient was used. For nanospray, a nanoElute UHPLC was fitted with an IonOpticks C18 column using a 210 minute gradient.

Results

The UPS1 standard was used to make a 5-step 1:3 dilution series in a constant background of NIST mAb over a concentration range from 0.3 to 934 ppm. PASEF enabled the detection of UPS1 proteins down to low single digit ppm concentrations in the presence of 25 µg NIST mAb. Linear response for the UPS1 protein beta-2-microglobulin in the concentration range from 132 ppm to 1.6 ppm was observed. The dilution series indicated that PASEF enabled detection of HCPs in the range of 1 to 100 ppm of the therapeutic protein. Nano UHPLC was also evaluated for its suitability for HCP identification. This setup has already been established as the new benchmark for bottom-up proteomics applications [1, 2]. More than 200 HCPs were identified in 1.5 µg NIST mAb, including expected and previously unreported HCP proteins. The depth of HCP coverage and sequence coverage was further extended to >280 HCP identifications by employing an alternate digestion method in which only the HCPs are digested allowing the mAb to be removed prior to analysis [3].

Conclusions

PASEF improved the sensitivity of peptide mapping enabling the detection of HCPs at sub 100 ppm levels. PASEF coupled to nanoLC facilitates detection of previously unreported trace level HCPs. The quality of MS/MS sequence spectra provided by the timsTOF Pro with PASEF allows high confidence in protein ID even when only 1-2 peptides are sequenced.

Novel Aspect

Sensitive HCP assay for routine use enabled by time and space focusing of precursor ions with a PASEF scan.

References

1. Meier et al., J Proteome Res., 14(12), 5378-87 (2015)
2. Lubeck et al., Bruker Application Note 131
3. Huang et al., Anal. Chem., 89, 5436-5444 (2017)

TP-50 / SIMULTANEOUS DETERMINATION OF TRAMADOL, O-DESMETHYLTRAMDOL, AND N-DESMETHYLTRAMADOL IN HUMAN PLASMA BY LIQUID CHROMATOGRAPHY-TANDEM MASS SPECTROMETRY AND ITS APPLICATION TO A PHARMACOKINETIC STUDY

Choong-Min Lee - Young-Hoon Kim - Won-Ki Chae - Eui Hyun Jung - Seok-Yong Lee

Sungkyunkwan University, School of Pharmacy, Suwon, Korea

Keywords: Tramadol; O-demethyltramadol; N-desmethyltramadol; HPLC-MS/MS; Human plasma; Pharmacokinetics

Introduction

Tramadol is a centrally acting analgesic drug and indicated for the treatment of moderate to moderately severe pain. To determine the tramadol and its two main metabolites, O-desmethyltramadol (ODT) and N-desmethyltramadol (NDT), in human plasma, a sensitivity-enhanced LC/MS/MS method was developed in this study.

Methods

A simple, rapid and sensitive high performance liquid chromatography analytical method coupled with tandem mass spectrometry (HPLC-MS/MS) for the simultaneous determination of the tramadol, ODT and NDT was developed and validated. Metoprolol was used as an internal standard (IS) and the analyte was extracted with methyl *tert*-butyl ether (MTBE). The chromatographic separation of tramadol, ODT, and NDT was achieved on a reversed-phase Luna C18 column (50 mm × 2.0 mm i.d., 5 μm particles), using a mixture of 10 mM ammonium formate buffer (pH 3.5)-methanol (30:70, v/v) at a flow rate of 200 μm/min. The total run-time was 2.5 min and the retention times of tramadol, ODT, NDT and IS were 0.88 min, 0.82 min, 0.90 min, and 0.85 min respectively. The analytes and IS were ionized using a positive ion mode and analyzed using multiple reaction monitoring(MRM) of the transitions at *m/z* 264.3→58.1 for tramadol, 250.3→58.2 for ODT, 250.3→44.1 for NDT, and 268.2→116.1 for IS.

Results

The weighted (1/x²) calibration curves were linear over plasma concentration range 0.20-400 ng/mL for tramadol (*r* = 0.999) and 0.25-200 ng/mL for ODT and NDT (*r* = 0.998 and *r* = 0.997, respectively). The lower limit of quantifications (LLOQ) of tramadol, ODT, and NDT using 100 μL of human plasma were 0.02 ng/mL, 0.25 ng/mL, and 0.25 ng/mL, respectively. This HPLC-MS/MS method showed improved sensitivity for quantification of tramadol, ODT, and NDT using lower volume of human plasma and reduced running time, compared with previously described analytical methods.

Conclusions

The method was successfully used to a pharmacokinetic study following the administration of a single 100 mg oral dose of tramadol to healthy male subjects.

Novel Aspect

To best our knowledge, there was only one case report using achiral determination of tramadol, ODT and NDT in human plasma employing HPLC-MS/MS and has lower sensitivity than this method

TP-51 / A RAPID AND HIGHLY SENSITIVE LIQUID CHROMATOGRAPHY-TANDEM MASS SPECTROMETRY METHOD FOR SIMULTANEOUS DETERMINATION OF MELOXICAM AND 5'-CARBOXYMELOXICAM IN HUMAN PLASMA AND ITS APPLICATION TO A PHARMACO

Choong-Min Lee⁽¹⁾ - ***Eui Hyun Jung***⁽¹⁾ - ***Yun Jeong Lee***⁽²⁾ - ***H. Y. Lee***⁽³⁾ - ***J. Lee***⁽¹⁾ - ***Seok-Yong Lee***⁽¹⁾

Sungkyunkwan University, School of Pharmacy, Suwon, Republic of Korea(1) - Dankook University, College of Pharmacy, Cheonan, Republic of Korea (2) - Konyang University, Dept. of Pharmacology, Daejeon, Republic of Korea (3)

Keywords: Meloxicam; 5'-carboxymeloxicam; LC-MS/MS; Human plasma; Pharmacokinetics.

Introduction

Meloxicam belongs to oxicam class of nonsteroidal anti-inflammatory drugs (NSAIDs) and is widely used to treat pain or swelling caused by osteoarthritis or rheumatoid arthritis. After oral administration of meloxicam, it is absorbed with bioavailability of 89% and has long half-life (22–24 hour). Meloxicam is extensively metabolized to four inactive metabolites by cytochrome P450 (CYP) 2C9 with a minor contribution of CYP3A4. The primary metabolic pathway of meloxicam is oxidation and 5'-carboxymeloxicam has been known for the major metabolite of meloxicam. The aim of the study was to develop and validate a more rapid and sensitive analytical method for simultaneous determination of meloxicam and 5'-carboxymeloxicam based on LC-MS/MS system.

Methods

After addition of 100 μL of 1 M hydrochloride acid, 200 μL aliquot of human plasma was extracted with an internal standard (piroxicam) by simple liquid-liquid extraction using methyl-*tert* butyl ether. The chromatographic separation of each analyte was achieved on a reversed phase Luna C18 column (50 mm \times 2.0 mm i.d., 5 μm particles) in a run time of 3 min, using a mixture of methanol and 10 mM ammonium formate buffer (pH 3.5) (80:20, v/v) at a flow rate of 200 $\mu\text{L min}^{-1}$. The tandem mass spectrometry was operated in multiple reaction monitoring acquisition mode using positive ion mode.

Results

This method was linear over the concentration range of 2–3000 ng/mL for meloxicam and 0.5–50 ng/mL for 5'-carboxymeloxicam, respectively. The lower limit of quantifications of meloxicam and 5'-carboxymeloxicam were 2 ng mL⁻¹ and 0.5 ng mL⁻¹, respectively.

Conclusions

The present analytical method showed better sensitivity and reduced running time than a previous study and successfully applied to a pharmacokinetic study after administration of 15 mg oral dose of meloxicam to healthy male subjects.

Novel Aspect

A rapid, simple, sensitive analytical HPLC-MS/MS method for the determination of meloxicam and 5'-carboxymeloxicam in human plasma was developed and successfully validated. Our study showed 4-fold more sensitivity for both analytes and reduced run time, compared to the only study for the simultaneous determination of meloxicam and 5'-carboxymeloxicam.

TP-52 / INVESTIGATION OF AMYLOID-BETA 1-42 REACTIVITY AND STABLE ADDUCT(S) FORMATION BY MASS SPECTROMETRY-BASED APPROACH

Marina Naldi⁽¹⁾ - *Angela De Simone*⁽²⁾ - *Daniele Tedesco*⁽¹⁾ - *Manuela Bartolini*⁽¹⁾ - *Vincenza Andrisano*⁽²⁾

University of Bologna, Department of Pharmacy and Biotechnology, Bologna, Italy (1) - University of Bologna, Department for Life Quality Studies, Rimini, Italy (2)

Keywords: LC-MS analysis; amyloid β 1-42 aggregation; CORM-3; Alzheimer's Disease

Introduction

Amyloid β -peptide 1-42 (A β 42) aggregation is a central event in the pathogenesis of Alzheimer's disease (AD). Carbon monoxide-releasing molecules (CORMs) are able to carry and release controlled amount of CO and are known to exert anti-inflammatory/apoptotic activities [1]. The study focus on the evaluation of the inhibitory properties of the carbon monoxide-releasing molecule 3 (CORM-3) towards A β 42 aggregation and the investigation of its mechanism of action.

Methods

A β 42 samples were prepared as previously reported [2]. CORM3/A β 42 adduct(s) formation and inhibition of A β 42 oligomerization were monitored by ESI-Quadrupole-Time of Flight (Q-ToF) flow injection analysis using reserpine as internal standard. Thioflavin T (ThT) assay [3] was used to monitor fibril formation, while circular dichroism (CD) spectroscopy was applied to follow A β 42 conformational changes [4].

Results

MS analyses showed that CORM3 can form stable adduct(s) with A β 42 monomers in a dose- and time-dependent manner. Confirming the evidence reported in a previous study [5], adduct formation involves histidine residues on A β 42 skeleton and the Ru^{II}(CO)₂ portion of CORM3. Inhibition of A β 42 aggregation strongly depends on CORM3/A β 42 ratio and was confirmed by both ThT assay and CD analysis. The latter offered further insights into the mechanism of action of CORM3, highlighting that CORM3 could prevent the conformational shift toward an amyloidogenic β -sheet rich conformer.

Conclusions

CORM3 emerged as a promising inhibitory agent toward A β 42 aggregation. The application of a MS-based method allowed confirming the formation of CORM3/A β 42 adduct(s) and highlighted CORM3 inhibitory activity. Combination of MS and spectroscopic data, such as fluorescence and CD, shed light on the mechanism of action of this compound. Furthermore, MS and ThT studies quantitatively contributed in identifying optimal CORM3/A β 42 ratio for best inhibitory effect.

Novel Aspect

The work highlighted CORM3 is endowed with a good inhibition capacity toward A β 42 aggregation, in addition to known anti-inflammatory properties.

References

1. Motterlini R., Clark J. E., Foresti R., C. J. Circ. Res., 90, 17, (2002).
2. Bartolini M., Bertucci C., Bolognesi M. L., ChemBiochem, 8, 2152 (2007).
3. LeVine H., Methods Enzymol., 309, 274 (1999).
4. Bartolini M., Naldi M., Fiori J., Anal. Biochem., 41, 215 (2011).
5. Santos-Silva T., Mukhopadhyay A., Seixas J. D., JACS, 133, 1192 (2011).

TP-53 / DOSE-AND TIME-DEPENDENT EFFECTS OF HEPATOTOXICANTS ON ADME PROTEINS IN 3D HUMAN LIVER SPHEROIDS BY PRM

Nathalie Selevsek⁽¹⁾ - *Witold E. Wolski*⁽¹⁾ - *Laura Kunz*⁽¹⁾ - *Henrik Cordes*⁽²⁾ - *Vanessa Baier*⁽²⁾ - *Lars Kuepfer*⁽²⁾ - *Olivia Clayton*⁽³⁾ - *Adrian Roth*⁽³⁾ - *Ralph Schlapbach*⁽¹⁾

Functional Genomics Center Zurich, ETH Zurich, Zurich, Switzerland (1) - Institute of Applied Microbiology, RWTH Aachen, Aachen, Germany (2) - Roche Pharma Research and Early Development, Roche Innovation Center Basel, Basel, Switzerland (3)

Keywords: ADME, 3D spheroids, drug toxicity, targeted proteomics, PRM

Introduction

ADME genes (ADME: absorption, distribution, metabolism, and excretion) encode for the involved phase I-II enzymes and transporters, and they respond adaptively to xenobiotic exposure. To predict in vivo drug response in ADME proteins we applied a Parallel Reaction Monitoring (PRM)-based targeted proteomics approach to quantify the proteins in 3D human liver spheroids challenged with different dose regimes of several hepato toxicants up to two weeks of treatment.

Methods

After isolation from drug-treated 3D liver spheroids, proteins were digested and analyzed by Liquid Chromatography - Mass Spectrometry (LC-MS) on a Q Exactive™ HF instrument (ThermoFisher Scientific) in PRM mode. Raw MS data were directly imported into Skyline for peptide/protein quantification. Since ADME proteins are structurally represented in PBPK models the measured changes in protein abundance may be directly integrated in the corresponding models.

Results

To design targeted PRM assays for the ADME proteins, we focused mainly on the ADME genes that are involved in the pharmacokinetics of the tested liver toxicants. MS coordinates such as precursor and retention time information were selected from spectral libraries built from MS/MS spectra acquired by data-dependent analysis experiments of protein digests isolated from pooled drug-treated 3D liver spheroids. In total, dozens of PRM assays were developed and optimized for the detection of the ADME proteins. Overall, we could identify and quantify with high confidence the different ADME proteins with either two or three proteotypic peptides across all the drug-treated 3D liver spheroids. We demonstrated high reproducibility in protein abundance between biological replicates across all time points and dose regimes. We showed that the abundance profile of ADME proteins were drug and dosage regime specific. The observed changes were integrated in the various PBPK models to describe cellular adaptation in response to xenobiotic exposure.

Conclusions

This allowed in particular to simulate the to-be-expected effect of multiple drug administration on drug pharmacokinetic over time. As such the presented approach presents an important tool to predict in vivo drug response in ADME proteins after multiple administration in a joint approach of model-based assay design and targeted proteomics.

Novel Aspect

We report quantitative and qualitative changes of ADME proteins in drug-treated 3D human liver spheroids by label-free PRM assays.

TPS-S02 - Single cell

Chairs: Peter Nemes, Giuliana Bianco

TP-54 / FUNCTION OF PIWIL3 IN MAMMALIAN OOCYTES

Minjie Tan⁽¹⁾ - **Helena van Tol**⁽¹⁾ - **Mirjam Damen**⁽²⁾ - **Tom Stout**⁽³⁾ - **Wei Wu**⁽²⁾ - **Bernard Roelen**⁽¹⁾

Utrecht University, Farm Animal Health, Utrecht (1) - Utrecht University, Biomolecular Mass Spectrometry & Proteomics, Utrecht (2) - Utrecht University, Equine Sciences, Utrecht (3)

Keywords: Affinity purification crosslinking mass spectrometry (AP-XL-MS), single-cell sequencing, PIWI

Introduction

PIWIs are argonaute proteins that interact with small non-coding RNAs (PIWI interacting RNAs; piRNAs). PIWI-piRNA complexes serve to target transposons for degradation, and are crucial guardians of genome integrity. PIWIL3 is specifically and highly expressed in unfertilized human oocytes, but a homologue is absent in mouse and rat, hence conventional animal systems cannot be used to study the critical function of PIWIL3 before fertilization[1].

Methods

We performed an in-depth analyses of bovine oocytes by mass spectrometry. From manually dissected oocytes, proteins were extracted by repeated freeze-thaw cycles and sonication. In parallel, we (i) performed AP-XL-MS on the bovine PIWIL3 complex, and (ii) profiled the oocyte proteome by LC-MS/MS on a Q-Exactive HF mass spectrometer, with pre-fractionation at high pH. We also followed the expression of PIWIL3 by microinjection and immuno-labeling techniques.

Results

The bovine system is the only large mammal system that expresses a PIWIL3 homologue. Herein, we investigated the expression, distribution, interaction, and argonaute function of PIWIL3 in bovine oocytes. By microinjection of PIWIL3-EGFP mRNA, we followed the localised expression of PIWIL3 in oocytes and early-stage fertilized embryos. From extremely limited oocyte material, we identified a total of 3897 proteins, which revealed the functional state of oocytes prior to fertilisation. In parallel, we identified a novel PIWIL3 interaction network involving two known regulators of argonaute function, which was independently verified by immuno-staining confirming intracellular co-localisation. Additional experiments involving single-cell sequencing and pi-RNA sequencing further support a critical role for mammalian PIWIL3 in protection of genome integrity prior to fertilization and generation of a new genome.

Conclusions

By a combination of AP-XL-MS and pi-RNA sequencing, we shed light on the role of PIWIL3 in mammalian oocytes prior to fertilization, and assign PIWIL3 to an interaction network in which its function may be further characterised.

Novel Aspect

We showcase here the extreme sensitivity of AP-XL-MS approach in mapping protein complexes from limited starting material.

References

1. Roovers et al., Cell Rep 2015;10, 2069-2082.

TP-55 / PICOLITER PRESSURE-PROBE-ELECTROSPRAY-IONIZATION MASS SPECTROMETRY FOR IN SITU MONITORING SINGLE PLANT CELL METABOLITE CHANGES UNDER STRESS CONDITIONS

Rosa Erra-Balsells⁽¹⁾ - ***Keisuke Nakata***⁽²⁾ - ***Hiroshi Wada***⁽³⁾ - ***Hiroshi Nonami***⁽²⁾

Departamento de Química Organica, Facultad de Ciencias Exactas y Naturales, Universidad de Buenos Aires, Buenos Aires (1) - Faculty of Agriculture, Ehime University, Matsuyama (2) - NARO Kyushu Okinawa Research Center, Fukuoka (3)

Keywords: pressure probe ESI, single cell, organelle, metabolomics.

Introduction

Cytoplasmic streaming occurs consuming the ATP energy formed during glycolysis, respiration and photosynthesis, moving organelles and metabolites all around plant cells. CPP is an instrument that can measure the water status and several properties of intact cells (*i.e.*, volume) [1]. The CPP-mass spectrometer combination (picoPPESI-MS) [2,3] allows to measure water status, cytoplasmic streaming velocities and metabolites in cells sequentially.

Methods

Picopressure Probe for water relations. Osmometer for cell turgor.

Internal Electrode Capillary picoPressure Probe Electrospray Ionization (picoPPESI) as ions source, with Orbitrap mass spectrometer (Exactive Plus, Thermo Fisher Scientific Inc., MA) as MS analyzer.

nanoESI- Orbitrap mass spectrometer (Exactive Plus, Thermo Fisher Scientific Inc., MA).

Results

Plants were subjected to water stress. Water status (turgor pressure; water potential) slightly decreased and they recovered by osmotic adjustment in 12 h. The cytoplasmic streaming velocity was maintained during osmotic adjustment. PicoPPESI-MS indicated that several phosphate-metabolites (m/z 78.96 [HPO_3^-], 96.97 [H_2PO_4^-], and 259.02 [Glu-6-phosphate]) were transiently accumulated in the cells in response to the treatment. Furthermore, not only metabolites yielded by glycolysis occurred in cytosol (m/z 87.01 [Pyruvic ac.], 179.06 [Hex], 259.02 [Glu-6-phosphate], and 341.11 [Hex2]) but also those from Krebs cycle in the mitochondrion (m/z 115.00 [Fumaric ac.], 117.02 [Succinic ac.], 133.01 [Malic ac.], 145.02 [α -Ketoglutaric ac.], 173.01 [Aconitic ac.] and 191.02 [Citric ac.]) accumulated in response to the stress treatment. Because these metabolic pathways are essential for ATP hydrolysis, ATP production and the osmotic adjustment, is obvious their contribution to maintain cytoplasmic streaming providing energy.

Conclusions

Tacking into account the metabolic pathways mentioned above and the results obtained by picoPPESI-MS monitoring the cell metabolites, ATP hydrolysis, ATP production and the osmotic cell adjustment are essential for keeping cytoplasmic streaming providing the necessary energy.

Novel Aspect

PicoPPESI mass spectrometry system has a potentiality to monitor metabolic changes in situ at cell level as well as at cell level.

References

1. Boyer J. S., Measuring the Water Status of Plants and Soils, Academic Press, San Diego (1995).
2. Gholipour Y., Erra-Balsells R., Nonami H., Anal. Biochem., 433, 70 (2013).
3. Nakashima T., Erra-Balsells R., Nonami H., Anal. Chem., 88, 3049 (2016).

TP-56 / DRUG UPTAKE AND METABOLISM IN SINGLE LIVER CELLS

Liliana Pedro⁽¹⁾ - **Patrick Rudewicz**

Novartis Institutes for BioMedical Research, GDC, Emeryville, United States

Keywords: single cell, drug metabolism, nanospray, high-resolution mass spectrometry

Introduction

Cells, even within a genetically identical culture, possess individual differences that cannot be detected by

analyzing cell populations. These individual differences can be crucial to cellular health and survival. Herein, we present a method that combines automated live single cell sampling with high resolution mass spectrometry that enables the detection and quantification of drug and drug metabolites in individual liver cells.

Methods

Amiodarone was used as a test compound. HepG2 and cryopreserved human hepatocytes were used as liver cellular models. Cells were incubated with 10 μ M amiodarone for 24h at 5% CO₂ and 37°C. Live cells were visualized in real-time during automated single cell sampling using a confocal microscope modified with a XYZ stage robot (Yokogawa Electric Corporation). Platinum-coated glass capillaries (cellomics tip, HUMANIX) were used for cell sampling and nanospray ionization. A Bruker static nanospray source and a 7T SolariX FT-ICR-MS were used for mass spectrometry analyses.

Results

CT-10 tips (10 μ m I. D.) were chosen to sample whole hepatocyte cells, which were maintained at 5% CO₂ and 37°C inside the sample vessel of the confocal microscope until sampled. Before automatic live cell sampling, series of slice images (z-stacks) were acquired to determine the volume of each cell. The cellular metabolism was quenched immediately after cell sampling. Prior to mass spectrometry analysis, an organic solvent mixture spiked with internal standard was loaded into the tip. Direct infusion of single cell contents resulted in complex mass spectra, with varied ion abundances and many interfering background ions. The resolving power of the SolariX FT-ICR-MS enabled unambiguous amiodarone identification and, operated in CASI mode, an improvement of at least 10-fold in sensitivity. The internal standard was used for both lock-mass calibration and intensity normalization. Concentrations of amiodarone and its major metabolite (*N*-desethylamiodarone) were integrated over total cell spray times and multiplied by the flow rates to afford the amount of each compound within each cell. Single cell concentrations were determined by normalizing these amounts to individual cell volumes. Single cell measurements agreed well with determinations of a hypothetical average cell obtained from cell lysates. Single cell measurements were, nonetheless, more variable than technical replicate analysis of cell lysates.

Conclusions

Pioneering work conducted by Masujima and colleagues had previously demonstrated the applicability of cellomics tips for both live single cell sampling and direct infusion into the mass spectrometer by nanospray ionization without prior sample processing [1,2]. In the current work, we confirmed the feasibility of the approach and addressed many technical challenges currently associated with single cell analysis: limited sample amount, high dynamism of cellular metabolic state, low sampling throughputs, high complexity of ion mixtures and varied abundances. With this method, absolute compound quantification in the actual cells was achieved.

Novel Aspect

A new method for determining compound uptake, accumulation and metabolic stability inside individual liver cells.

References

1. Mizuno H., Tsuyama N., Harada T., Masujima T., *Journal of Mass Spectrometry*, 43, 1692-1700 (2008).
2. Fujii T., Matsuda S., Tejedor M. L., Esaki T., Sakane I., Mizuno H., Tsuyama N., Masujima T., *Nature Protocols*, 10, 1445–1456 (2015).

TP-57 / THE ASSEMBLY OF ISLET AMYLOID POLYPEPTIDE (IAPP) AND THE ROLE OF THE DISULFIDE BOND

Michael Bowers - Alexandre Ilitchev - Dezmond Bishop - Steven Buratto

UCSB, Chemistry and Biochemistry, Santa Barbara

Keywords: Ion Mobility, Atomic Force Microscopy, IAPP, Amyloid

Introduction

Islet Amyloid Polypeptide (IAPP), a 37 residue hormone, has been identified as a primary agent in the etiology of Type 2 Diabetes and the driving force for pancreatic cell death. There has been a substantial amount of research on the mechanism of the cytotoxic process but it is still not well understood. It is known that assembly of IAPP is involved and the oligomer states not the fibrils are the toxic agents. Here we provide new evidence bearing on the molecular details of the IAPP assembly process.

Methods

Ion Mobility Based Mass Spectrometry (IMS-MS) will be used to characterize the assembly of IAPP and a number of N-terminal mutants designed to shed light on the importance of the disulfide bond in this process. These results will be complimented with data from Atomic Force Microscopy (AFM) using the same solutions utilized in the IMS-MS studies.

Results

There is a disulfide bond linking cysteine residues at positions 2 and 7 in wild type IAPP. The stability of this loop appears to limit the participation of the N-terminal regions in the amyloid assembly of IAPP. However recent studies have shown that the N-terminal region is capable of self-assembly. Further, wt IAPP has an inhibitory affect on insulin secretion in beta cells but truncating the N-terminus reduces this effect. Here we look at a series of modified IAPP molecules and compare their assembly with wt: These include reduced IAPP, IAPP(8-37, and cysteine side chain Carbamidomethylation (CAM). We find reduction of the disulfide bond results in significant acceleration of the aggregation process while the CAM mutant greatly inhibits aggregation. The truncated form behaves similar to wild type in formation of oligomers but accelerates fibril formation. The implication of these results will be discussed.

Conclusions

The N-terminus has been shown to have a dramatic effect on the assembly of IAPP contrary to expectations. The results presented here contribute to clarifying the assembly mechanism of IAPP which is crucial to developing therapeutic agents to ameliorate its cytotoxic effect in Type 2 Diabetes

Novel Aspect

New data on the importance of the N-terminal region of IAPP is presented and the complimentary nature of IMS-MS and AFM emphasized.

TP-58 / GAS-PHASE COMPLEXATION OF A-/B-CYCLODEXTRIN WITH AMINO ACIDS STUDIED BY ION MOBILITY-MASS SPECTROMETRY AND MOLECULAR DYNAMICS SIMULATIONS

Yinjuan Chen - Zhicheng Zuo - Wenning Wang - Chuanfan Ding

Fudan University, Department of Chemistry, Shanghai

Keywords: Non-covalent Complexes; cyclodextrin; amino acids; ion mobility-mass spectrometry; molecular dynamics simulation

Introduction

Cyclodextrins (CDs) have the hydrophilic surface and hydrophobic cavity and enable CDs to form non-

covalent complexes with other molecules. Electrospray ionization-mass spectrometry (ESI-MS), and ion mobility-mass spectrometry (IMMS) are suitable technology for the gas phase conformational study. Herein, non-covalent complexes involving amino acids and α -/ β -CD are investigated using IMMS combined with molecular dynamics simulations.

Methods

The gas-phase complexation of non-covalent complexes involving α -/ β -cyclodextrin (α -/ β -CD) and amino acids (AA) (AA=Gly, Leu and Phe) are investigated by ion mobility-mass spectrometry in negative mode. The complex ions are further identified by tandem mass spectrometry. The complex gas-phase conformations are clarified by IMS in combination with molecular dynamics simulations.

Results

The 1:1 [α -/ β -CD+AA]⁻ adducts are observed when equimolar α -/ β -CD:Gly:Leu:Phe solutions are electrosprayed. The complex ions are further identified by tandem mass spectrometry, in which complex ions are fragmented into [CD]⁻ via losing one neutral AA molecule. Their gas-phase stability order of [α -/ β -CD+Gly]⁻ < [α -/ β -CD+Leu]⁻ < [α -/ β -CD+ Phe]⁻, revealed by survival yield curves, can be explained by gas-phase basicity difference. The anion is kinetically more stable with much smaller gas-phase basicity difference of [CD]⁻ and [AA]⁻. The complex gas-phase conformations are clarified by IMS in combination with molecular dynamics simulations. The nearly equal collision cross sections of [β -CD+Gly/Leu/Phe]⁻ anions were obtained by ion mobility analysis, while various cross sections were observed for the [α -CD+Gly/Leu/Phe]⁻ anions.

Conclusions

The structure of [α -CD+Leu/phe]⁻ which obtained by IMS represented the exclusion structure in which the hydrophobic side chain remains outside the hydrophobic cavity of α -CD, and inclusion structure was proposed for other complex anions with the AA molecule being trapped inside the hydrophobic cavity of CD. The conformations are consistent with the most stable configurations obtained by molecular dynamics simulations.

Novel Aspect

The gas phase inclusion or exclusion structure of non-covalent complexes of α -/ β -cyclodextrin and amino acids can be identified by IMS and molecular dynamics simulations.

References

1. P. Xiao, H. Li, M. Li, D. Song, X. Li, X. Dai, Z. Hu, Anal. methods 7 (2015) 8432-8438.
2. W. Wei, Y. Chu, R. Wang, X. He, C. Ding, Rapid Commun. Mass Spectrom. 29 (2015) 927-936.
3. M. Rejzek, C. E. Stevenson, A. M. Southard, D. Stanley, K. Denyer, A. M. Smith, M. J. Naldrett, D. M. Lawson, R. A. Field, Mol. Biosyst. 7 (2011) 718-730.
4. G. Polekhina, A. Gupta, B. van Denderen, S. C. Fell, B. E. Kemp, D. Stapleton, M. W. Parker, Structure 13 (2005) 1453-1462.

TP-59 / INVESTIGATING THE INTERLOCKED NATURE OF CATENANES AND KNOTS WITH MASS SPECTROMETRY AND ION MOBILITY SPECTROMETRY

Anneli Kruve-Viil - Christoph A. Schalley

Freie Univeversität Berlin, Institut für Chemie und Biochemie, Berlin

Keywords: mechanical bond, ion mobility, catenanes, knots

Introduction

Mechanically interlocked molecular architectures are exclusively used to convert energy into directed molecular motion on molecular level. Catenanes [1-3] and other mechanically interlocked systems are often hard to unambiguously characterize; the NMR studies need isolation and purification of the MIM as well as control experiments. On the other hand, in mass spectrometry (MS) unspecific binding often cannot be ruled out.

Methods

We are proposing a fast, selective, and simple approach based on collision induced dissociation combined with ion mobility spectroscopy (CID/IMS) to unambiguously characterize the topology for catenanes and knots.

Results

CID is used to induce the cleavage of the covalent bonds and IMS to separate the formed species based on their size in gas phase. Therefore, applying CID on catenanes results on two species an open macrocycle and a closed macrocycle, which can be easily separated with IMS and detected with mass spectrometry. The method automatically involves a quantitative control experiment: the amount of formed open and closed macrocycle sum up to the amount of catenane that has under gone cleavage. Similarly, this technique helps distinguishing knots from simple macrocycles. We have tested CID/IMS based method to characterize several libraries with already well known architectures including imine-bond based catenanes in water, oligolactame based structures, metal templated catenanes, etc. Thereafter, we applied CID/IMS method to probe the libraries with previously unidentified topology. Both catenanes and knots were identified; however, also libraries without mechanically interlocked bonds could be identified as true negatives. For the unknown libraries, the structures were later separated and confirmed also with other techniques.

Conclusions

We demonstrate that CID/IMS based method is rapid, conclusive, straightforward, and can be applicable to various types of catenane, including both for symmetrical and unsymmetrical catenanes, mixtures, or even crude products. Additionally, in case of dynamic libraries it can also be used to effectively probe the kinetic and thermodynamic properties of the mechanical bond.

Novel Aspect

MS, combined with IMS and IMS/MS/MS can be effectively used to characterize the structure of interlocked supramolecular structures.

References

1. Caprice K., Pupier M., Krueve A., Schalley C. A., Coughon F., *Chemical Science* 9, 1317 (2018).
2. Schalley C. A., Hoernschemeyer J., Li X., Silva G., Weis P., *International Journal of Mass Spectrometry* 228, 373 (2003).
3. Browne C., Ronson T. K., Nitschke J. R., *Angewandte Chemie Internationala Edition* 53, 10701 (2014).

TP-60 / COLLISION INDUCED DISSOCIATION OF MULTIPLY CHARGED PEPTIDE AGGREGATES

Mohammed Al-Majidi⁽¹⁾ - ***Anna Magyar***⁽²⁾ - ***Gitta Schlosser***⁽³⁾

ELTE Eötvös Loránd University, Doctoral School of Chemistry, Budapest (1) - Hungarian Academy of Sciences, ELTE Eötvös Loránd University, MTA-ELTE Research Group of Peptide Chemistry, Budapest (2) - ELTE Eötvös Loránd University, Department of Analytical Chemistry, Budapest (3)

Keywords: clusters, charge separation, non-covalent complexes, peptides, tandem mass spectrometry

Introduction

Self-association and interactions of amino acids and peptides is in the focus of intensive research since a long time. In this work, clustering properties and gas-phase dissociation pathways of multiply charged peptide clusters were studied. We have analyzed the cluster formation of tripeptides as well as the collision induced dissociation of their clusters in an ion trap mass spectrometer.

Methods

Mass spectrometric experiments were performed using two tyrosine-containing tripeptides: H-Tyr-Ala-Gly-OH and H-Tyr-Gly-Phe-OH, on a Bruker Esquire 3000+ ion trap mass spectrometer equipped with ESI source. Optimized ion source parameters were: nebulizer gas pressure: 80 psi, dry gas flow: 0.5 l/min, dry temperature: 250°C, skimmer voltage: 40V, capillary exit voltage: 106 V.

Results

The relatively high peptide concentrations were found to enhance the formation of peptide aggregates. Ion source parameters were optimized to achieve high signal intensity for larger peptide aggregates. We observed, that high nebulizer gas pressure and low dry gas flow rate facilitate the detection of large, multiply charged clusters. Both peptides exhibit a characteristic cluster distribution profile in ESI-MS. Beside the protonated monomer, dimer and trimer species, peptides produce larger clusters in the 1200-2400 m/z . Most intensive peaks correspond to 5 and 6 peptide units with one charge. MS/MS experiments revealed that multiply charged metaclusters of peptides appear at the same mass-to-charge ratios as peptide aggregates with lower charge states. Peptide clusters of interest were isolated in the ion trap and fragmented to identify their fragmentation pathways.

Conclusions

The dominant fission products for the tripeptide clusters were the losses of the proton-bound dimers, presumably due to the high gas-phase stability of this species. Charge separation by cluster fission was proved to be a major dissociation pathway of these peptide clusters. Beside the formation of singly charged ion population, triply charged species were preferred products as well.

Novel Aspect

Characteristic cluster distribution was detected for tripeptides, together with typical charge separation-based fragmentation pathways.

Acknowledgement

Gitta Schlosser acknowledges the support of the MTA Premium Post-Doctorate Research Program of the Hungarian Academy of Sciences (HAS, MTA).

TP-61 / NOVEL, BILATERAL FLUOROUS ALKOXIDE-BASED CLUSTERS DISCOVERED BY AMBIENT IONIZATION MASS SPECTROMETRY

Máté Berta⁽¹⁾ - **Tibor Nagy**⁽²⁾ - **Sándor Kéki**⁽²⁾ - **József Rábai**⁽¹⁾ - **Gitta Schlosser**⁽³⁾

Department of Organic Chemistry, ELTE Eötvös Loránd University, Budapest (1) - Department of Applied Chemistry, University of Debrecen, Debrecen (2) - Department of Analytical Chemistry, ELTE Eötvös Loránd University, Budapest (3)

Keywords: fluorous flying salts, fluorous clusters, ambient ionization, non-covalent complexes, energy-dependent tandem mass spectrometry

Introduction

Cluster formation is an important research topic in fundamental mass spectrometry. Formation of alkali

metal halide clusters using electrospray ionization has long been known. These salts can form very large cluster structures such as $\text{Na}_{124}\text{Cl}_{122}^{2+}$ [1] Here we report the discovery of novel, bilateral fluorine-based clusters by ESI, showing unique structural properties.

Methods

Experiments were performed using a Bruker Esquire 3000+ ion trap MS with ESI ion source and a Bruker microQTOF with DART ion source. MM/QM computational studies were carried out with Schrödinger Maestro 11.1 and Gaussian16 software packages. We optimized the structures with OPLS2005 force field and B3LYP/Def2SVP (GD3BJ) and calculated single point energy with B3LYP/Def2TZVPP level of theory.

Results

Nonafluoro-*tert*-butyl alkoxides, $\text{MOC}(\text{CF}_3)_3$ [$\text{M} = \text{Li}, \text{Na}, \text{K}, \text{Rb}, \text{Cs}$], were synthesized to investigate their clusters using mass spectrometry. Solutions prepared in dry acetonitrile were directly injected into the ESI ion source. The salts show complexes of outstanding stability and high intensities both in positive and negative ionization modes, such as $[\text{M}_4\text{A}_3^+; \text{M}_5\text{A}_6^-; \text{where } \text{A} = (\text{OC}(\text{CF}_3)_3)^-]$.

Computational chemistry methods, MM/QM modelling, were used to calculate the stable gas phase structures. Energy-dependent MS/MS fragmentation experiments were also performed to quantitatively describe affinity and stability parameters of clusters. Experimental data were in good correlation with calculated values.

Conclusions

Formation of fluorine-based alkoxide clusters are described here for the first time. These special compounds, referred to “fluorine flying salts”, represents a new cluster type, in which the anion has fluorine-ionic bilateral properties, which prevents the limitless polymer cluster formation, and support the outstanding stability of defined structures.

Novel Aspect

Bilateral fluorine-based alkoxide clusters discovered by ESI-MS give novel opportunity for the study of electron pair acceptor compounds by MS.

References

1. Hao C., et. al, J. Mass Spectrom., 36, 79-96 (2001).

Acknowledgement

Gitta Schlosser acknowledges the support of the MTA Premium Post-Doctorate Research Program of the Hungarian Academy of Sciences (HAS, MTA).

TP-62 / NONCOVALENT COMPLEXES IN NATURAL EXTRACTS

Pietro Traldi⁽¹⁾ - **Sara D'Aronco**⁽²⁾ - **Sara Crotti**⁽¹⁾ - **Marco Agostini**⁽²⁾ - **Luisa Mattoli**⁽³⁾ - **Michela Burico**⁽³⁾ - **Sara Tamimi**⁽³⁾

Fondazione Istituto di Ricerca Pediatrica Città della Speranza, Città della Speranza, Padova (1) - Dipartimento di Scienze Chirurgiche, Oncologiche e Gastroenterologiche, Università di Padova, Padova (2) - Aboca, Ricerca Fitochimica, Sansepolcro (3)

Keywords: noncovalent complexes in natural extracts

Introduction

Chemical evolution includes capture, mutation and propagation of molecular information via covalent and noncovalent associations [1]. Considering the high complexity of natural substrates, for the presence of

molecular species of different chemical nature, the formation of noncovalent complexes has to be expected. They could play a role in the activity of the extract, leading to structures of active compounds different to that of the synthetic analogous products.

Methods

The presence of noncovalent complexes in different natural extracts (green tea, Ceylon tea, green coffee, mate) has been investigated by precursor ion scan on ions generated by ESI and neutral loss scans. In order to exclude that the complexes so detected were artifact due to ESI conditions, LC MS measurement were performed, and reconstructed chromatograms based on accurate mass values confirmed the presence of noncovalent complexes in the natural extracts.

Results

Some investigations were performed to establish unequivocally the presence of catechin/caffeine complexes in green tea (GTE), Ceylon tea, green coffee, mate extracts. By ESI-MS experiments performed on GTE by means of precursor ion scan and neutral loss scan methods, the presence of caffeine:catechin complexes was proved [2]. The same approach showed that in the case of Ceylon tea the abundance of these complexes is strongly reduced, while for green coffee and mate complexes of caffeine with chlorogenic acid present in the extracts were undetectable. Further experiments were performed on GTE by LC/MS method operating in high resolution conditions. The reconstructed ion chromatograms of the exact mass ions corresponding to the caffeine/catechin complexes have been obtained, showing the presence of the complexes of caffeine with gallate-type catechins [2]. This approach was able to put in evidence the presence of the same complex with different structures, exhibiting different retention times, in agreement with the data of Ujihara and Hayashi [3].

Conclusions

The complexes between caffeine and catechins already observed in solid state and in solution by NMR in ad hoc prepared solution are detectable also by mass spectrometry, once specific methods are employed. These complexes are detected in green tea and, in minor extent, in Ceylon tea. In mate and in green coffee the complexes possibly originating by interaction of caffeine with chlorogenic acid are undetectable, reasonably due to the low affinity of the two molecular species.

Novel Aspect

For the first time the presence of caffeine/catechin complexes in green tea and in Ceylon tea extracts has been proved by mass spectrometric experiments.

References

1. Linn D., Burrows C., Goodwin J., Metha A., Acc. Chem. Res. 45, 2023-2024.
2. Mattoli L., Mercati V., Burico M., Porchia M., Tisato F., D'Aronco S., Crotti S., Agostini M., Traldi P., J. Nat. Prod., submitted
3. Ujihara, T. ; Hayashi, N. J. Nat. Prod., 2016, 79, 66-73

TP-63 / STERIC STRUCTURE AND STABILITY OF THE GAS-PHASE COMPLEXES OF MACROCYCLIC RECEPTORS WITH ANIONS

Magdalena Zimnicka - Alicja Wasilewska - Kajetan Dabrowa - Witold Danikiewicz

Polish Academy of Sciences, Institute of Organic Chemistry, Warsaw, Poland

Keywords: ion mobility mass spectrometry, macrocycles, anions, noncovalent complex

Introduction

One of the most intensively developed topics of supramolecular chemistry is design, synthesis and studies

of neutral molecular receptors capable of strong and selective binding of anions because of their crucial role in many chemical and biochemical processes. The ion mobility mass spectrometry (IM-MS) supported by theoretical calculations and CID methods are evaluated as fast screening methods in the development of new receptors.

Methods

IM-MS measurements were performed on Synapt G2-S HDMS. Collision-induced dissociation (CID) measurements were performed on the triple quadrupole 4000 Q Trap (MD SCIEX). Starting models of the molecules and ions were generated using Spartan software and were further optimized with DFT methods. The theoretical collision cross sections for model structures were calculated using Mobcal program.

Results

The complexes of macrocyclic receptors with selected anions are stable enough in the gas-phase to survive the conditions of the IM-MS experiments. The experimentally measured collisional cross sections were used for establishing the steric structures of the studied ions by comparing them with the cross section values obtained from molecular modelling methods. The type of calibrant ions significantly affects the accordance between theoretical and experimental values of collision cross sections. The effects of the deformation of macrocyclic receptors upon complexation with anions are recognized in IM-MS experiments.

Conclusions

The results confirm the usefulness of the IM-MS technique for studying the structures of complexes of the macrocyclic receptors with anions. They also prove that the IM-MS method is very sensitive to even small differences in the structures of the studied complex ions. The IM-MS and CID techniques are potentially very efficient and cost-effective methods to study supramolecular recognition of anions.

Novel Aspect

IM-MS and CID techniques evaluated as fast screening methods in the development of new receptors.

TPS-S04 - High resolution MS: instrumentation, methods & applications

Chairs: Evgeny Nikolaev, Gianluca Giorgi

TP-64 / TOWARDS ROUTINE INTACT LEVEL QUANTIFICATION OF MONOCLONAL ANTIBODY IN PLASMA BY HRMS

Yun Alelyunas - Mark Wrona - Henry Shion - Caitlin Dunning - Jayne Kirk

Waters Corp, Scientific Operations, Milford

Keywords: HRMS, Monoclonal Antibody, Quantification, Intact, Bioanalysis

Introduction

Monoclonal antibodies (mAbs) are an important class of innovative medicine for the treatment of cancer and other immune mediated diseases. Five of the top ten best-selling medicines in 2017 worldwide were mAbs. Traditionally, bioanalysis of mAbs has been done almost exclusively using ligand binding assays with some structural characterization provided by high resolution mass spectrometry (HRMS). More recently, direct intact level quantification based on HRMS has gained attraction, as it offers direct analysis, simpler sample preparation compared to surrogate peptide approach, and also provides some structural and isoform information associated with the mAb in study. In this presentation, the quantification of several

mAbs in biological matrix is performed. Best practices for sample preparation, LC separation, MS detection, and data processing for achieving the high sensitivity are described. Challenges and the needs for next generation sample preparation and analysis of intact mAb analysis are discussed.

Methods

Adalimumab, infliximab, or trastuzumab samples were prepared by spiking mAb solution into mouse plasma. Each mAb was immuno purified from plasma using goat anti-human Fc prepared either in streptavidin coated plate (1) or magnetic beads. A Tecan Genesis liquid handling platform was employed in sample preparation. After elution under acidic conditions, extracted mAb was quantified by UPLC-HRMS. The LC-MS system consists of a binary pump, sample manager, column manager, and a Vion QToF HRMS. The mass range used for data acquisition was 500-4000 m/z , and the system was operated in ESI+ full scan sensitivity mode. A generic gradient was used with a flow rate of 0.2 mL/min and a total run time of 13 minutes. Data was acquired, analyzed and quantified using the UNIFI informatics system.

Results

A standard calibration curve was prepared by spiking a serial dilution of monoclonal antibody into mouse plasma using BSA as a carrier protein. After immunoaffinity purification, the solutions were analyzed by LC-MS under generic reversed phase conditions using LC/MS friendly mobile phase consisted of acetonitrile, H₂O, and formic acid. LC-MS conditions were optimized including chromatographic separation of the mAb from the carrier protein, glycan isoform spectral profile, and charge state summation and XIC tolerances for processing, yielding best bioanalytical performance. Excellent sensitivity (<50 ng/mL) using the optimized conditions was achieved. Comparing immunoaffinity purification from streptavidin plate vs magnetic bead, it was found that the plate offered lower background signal compared to magnetic bead. Preliminary results showed that LLOQ determined from streptavidin coated plate is <50 ng/mL compared to 250 ng/mL using streptavidin coated magnetic bead. Data processing parameters were also extensively evaluated and optimized. It was found that both XIC mass tolerance window and the number of charge states included in the peak area summation can have a large effect on data quality, especially at the LLOQ level. A flexible informatics workflow evaluating multiple parameters was used for a quickly viewing the data and determining the best analysis method.

Conclusions

HRMS quantification of intact mAbs offers analytical scientists many advantages over traditional approaches. In addition to simpler sample preparation (versus surrogate approaches), structural information at the whole molecule level can provide additional information in understanding the state of mAb in question in biological systems and interpreting quantitative outcomes. In this study, high sensitivity quantification of monoclonal antibody is presented. results show that the combination of a sensitive mass spectrometer with the best practice in sample preparation and data processing parameters are critical in obtaining the highest sensitivity. Combination of these best practices brings us a step closer to routine intact level quantification of mAbs in plasma.

Novel Aspect

The present study provides a complete workflow for intact level mAb quantification in plasma from automated sample preparation through automated data processing and data review. The LLOQ of <50 ng/mL represented the best sensitivity so far and further closing the gap between intact level and surrogate peptide approaches in the quantification of monoclonal antibody in plasma.

References

1. Kellie J., Kehler J., Mencken T. J., Snell R. J., Hottenstein C. S. *Bioanalysis* 8, 2103 (2016).

TP-65 / PROGRESS TOWARDS IMPLEMENTING ACCURATE MASS MS FOR ROUTINE BIOTHERAPEUTIC ANALYSIS

Henry Shion⁽¹⁾ - William Alley⁽¹⁾ - Jonathon Pugh⁽²⁾ - Robert Lewis⁽²⁾ - Ying Qing Yu⁽¹⁾ - Weibin Chen⁽¹⁾

Waters Corp, Scientific Operations, Milford (1) - Waters Corp, Waters, Wilmslow (2)

Keywords: fit-for-purpose time-of-flight MS, biotherapeutics, intact mass analysis, peptide mapping, glycan analysis

Introduction

The adoption of accurate mass MS for routine analysis has been hindered by instrument size, operation complexity, hardware robustness, and reproducibility. Yet, there is a need for this approach; for example, efficient data collection and analysis of quality-relevant information within the pharmaceutical industry. Here, we describe the performance evaluation of a prototype small-footprint time-of-flight MS that could overcome many of the above challenges.

Methods

To assess the performance and reproducibility across a range of biotherapeutics, trastuzumab, NIST mAb, infliximab, trastuzumab emtansine, brentuximab vedotin, and EPOs were analyzed using the prototype MS coupled with various chromatography separation modes for intact, subunit, peptide mapping and released glycan experiments. Data acquisition, processing, and reporting were all achieved within a single workflow-centric informatics platform.

Results

For intact mass analysis, a dilution series for both trastuzumab and NIST mAb established a limit of detection (LOD) and intra-system LOD reproducibility. Results were comparable to other high resolution MS systems. More structurally complex proteins (infliximab, trastuzumab emtansine, and EPOs), were analyzed and data demonstrated the greater heterogeneity of these molecules. The average drug-to-antibody ratio (DAR) for a cysteine-conjugated ADC agreed with values from commercial MS systems and hydrophobic interaction chromatography (HIC).

NIST mAb subunits showed good mass accuracy. Low oxidation levels were observed for both subunits after forced degradation experiments.

Peptide mapping of the NIST mAb was used to evaluate system robustness. A high sequence coverage was determined through precursor accurate masses and identification of fragment ions.

Glycan analysis was performed on the NIST mAb. All major glycans were correctly identified as well as many lower abundance glycans, including some potentially immunogenic glycans.

Conclusions

Several key biopharmaceutical experiments were performed on a small-footprint prototype TOF MS instrument. The data generated for all experiments (intact mAb mass, peptide mapping, glycan analysis) were very comparable to currently available commercial high resolution MS instrumentation. Thus, these data show great promise that this instrument could be an appropriate choice for routine MS-based biopharmaceutical analyses.

Novel Aspect

A prototype small-footprint oa-TOF MS provides fit-for-purpose data for multiple experiment types relevant to routine biotherapeutic analysis.

TP-66 / MULTIPLEXED ANALYSIS OF ENVIRONMENTAL SAMPLES BY GCXGC-HRMS WITH MULTI-IONIZATION CAPABILITIES

Viatcheslav Artaev⁽¹⁾ - *Georgy Tikhonov*⁽¹⁾ - *Albert Lebedev*⁽²⁾

LECO Corporation, R&D, St. Joseph (1) - Moscow State University, Chemistry Dept., Moscow (2)

Keywords: GCxGC, TOFMS, HRMS, Non-targeted, Environmental

Introduction

The Gas Chromatography coupled to Mass Spectrometry provides powerful analytical tools, combining selectivity, sensitivity, reliability, and information capacity for targeted or non-targeted methods of environmental analysis. Yet, the complexity of the environmental samples and variety of the analytes on interest in them demands the novel methods of analysis, taking advantage of modern advances in gas chromatography, mass spectrometry and data processing.

Methods

The GCxGC equipped with thermal modulator and coupled to High Resolution TOFMS was used for untargeted screening of snow samples. The mass spectrometer was equipped with a novel ion source, capable of three modes of ionization - EI, PCI and ENCI. The snow samples were analyzed in different ionization modes using same instrument and at the same chromatographic conditions. The chromatograms were aligned for assisting in identification of the analytes.

Results

The high resolution TOFMS is an ideal GC-MS detector for screening and quantifying unknowns as it collects high mass accuracy full mass range mass spectra at very high acquisition rates without any data loss. The technique provides highly reliable data suitable for automatic accurate spectral deconvolution of the coeluting analytes present in the samples in the wide concentration range. The GCxGC technique increases separation power, allows chromatographic separation of closely eluting constituents, making analytes identification more reliable and makes comprehensive analysis of environmental samples more realistic. Comparison of analytes detected in the same sample using EI, PCI and ENCI ionization shows that individual ionization method may miss certain types of the analytes, and combining data from three modes of ionization improves the quality of the results. The presentation will demonstrate the practical implementation of the technique and will report examples of false negative results if only one ionization method is used.

Conclusions

Application of the highly efficient multiplexed technique of GCxGC-HRTOF MS with multi-ionization capabilities has proved to be very powerful for increasing coverage and reliability of the untargeted screening results of environmental samples. Combining data from all three ionization modes increased the number of detected analytes by at least 30% and improved reliability of the analytes assignments.

Novel Aspect

GCxGC-HRMS analysis with EI, PCI and ENCI ionization of analytes was used for comprehensive analysis of snow samples collected at the same day in various parts of large city.

TP-67 / GAS CHROMATOGRAPHY-ATMOSPHERIC PRESSURE PHOTOIONIZATION-HIGH RESOLUTION MASS SPECTROMETRY FOR THE DETERMINATION OF SEMI-VOLATILE FLUORINATED COMPOUNDS

Juan Francisco Ayala-Cabrera - Encarnación Moyano - Francisco Javier Santos Vicente

University of Barcelona, Department of Chemical Engineering and Analytical Chemistry, Barcelona

Keywords: gas chromatography, atmospheric pressure photoionization, high resolution mass spectrometry, semi-volatile fluorinated compounds, water

Introduction

Semi-volatile fluorinated compounds are currently determined by GC-MS using electron and chemical ionizations [1], but problems related to their low ionization efficiency have been reported. Considering the proved advantages of atmospheric pressure photoionization (APPI) for the analysis of non-polar

compounds [2], in this work we explored the applicability of this technique for the determination of several families of fluorinated compounds by GC-MS.

Methods

AGC-APPI-HMRS (Orbitrap) method combined with HS-SPME was developed for the determination of fluorotelomer olefins (FTOs) fluorotelomer alcohols (FTOHs) and fluorinated sulfonamides (FOSAs) and sulfonamide-ethanols (FOSEs) in river water samples at trace levels. The HS-SPME method was performed using a PDMS/DVB/CARB fiber and the GC-APPI-HRMS determination was carried out in the negative ion mode using full scan acquisition at a resolution of 35,000 (FWHM).

Results

Semi-volatile fluorinated compounds were ionized via dopant-assisted in negative ion mode, generating the $[M+O_2]^-$ as base peak of mass spectra for most of FTOHs and FOSEs, while the $[M-H]^-$ and several in-source CID fragment ions were mainly observed for FOSAs and FTOs, respectively. Acetone was found to be the best dopant for achieving the highest responses for all the compounds. The GC-APPI-HRMS method provided low limits of detection, ranging from 0.05 to 10 pg injected, and non-significant matrix effects for the analysis of water samples with different pollutant loads. HS-SPME was used to extract analytes from water river samples and the method limits of detection achieved ranged from 24 $\mu\text{g L}^{-1}$ to 15 ng L^{-1} . The HS-SPME GC-APPI-HRMS method provided good linearity ($R^2 > 0.993$) in the calibration working range. The precision and trueness of the developed method were determined using spiked blank river water samples, obtaining for all the compounds RSD% values lower than 11% and relative errors down to 12%, respectively.

Conclusions

GC-APPI-HRMS has shown to be an excellent alternative to the existing GC-MS methods using classical ionizations for the determination of semi-volatile fluorinated compounds, providing LODs up to 80 times lower than EI and CI techniques. Moreover, the response of the GC-APPI-HRMS is not significantly affected by the matrix. On that way, the developed GC-APPI-HRMS method in combination with HS-SPME allowed the detection of the target compounds down to $\mu\text{g L}^{-1}$ levels.

Novel Aspect

This is the first time that a GC-APPI-HRMS method is proposed for the determination of semi-volatile fluorinated compounds in river water samples at very low concentrations.

References

1. Szostek B., Prickett K. B., Buck R. C., Rapid. Commun. Mass Spectrom., 20, 2837 (2006).
2. Li D-X., Gan L., Bronja A., Schmitz O. J., Anal. Chim. Acta, 891, 43 (2015).

TP-68 / DETERMINATION OF QUANTITATIVE AMOUNTS AND RATIO BETWEEN ENZYME VARIANTS IN A SAMPLE, BY SPIKING WITH STABLE ISOTOPE LABELLED PEPTIDES

Janne Bech Thøgersen

DuPont Nutrition Biosciences ApS, Advanced Analysis, 8220 Brabrand

Keywords: Protein mutation variants, quantification, isotopic labelled peptide standards

Introduction

A sample was found to contain two similar protein variants ("F" and "Y"), only differing by amino acid substitutions at a few positions. An attempt to estimate the ratio between the two variants had previously been carried out, using proteolytic digestion of the sample. Two peptides separating the variants were used

for calculation of a chromatographic peak-area-ratio. However, this method suffered from poor precision mainly caused by variations in the nanoLC-MS system. Accordingly, the aim of the present work was to improve the precision of the ratio determination and to know the exact quantitative amount of the variants in the sample.

Methods

Initially, protein concentration was estimated using a BCA assay, to control protein to enzyme ratio, and for the final quantitative calculations of the content of variants relative to the total protein amount.

One of the positions at which the variants were differing, was captured by digestion using endo-proteinase Asp-N, giving two good responding peptides. The Asp-N digested sample was spiked with corresponding AQUA peptides, purchased as synthesized variants with the C-terminal amino acid being heavy isotopically labelled (+4 Da). Accordingly, it was possible to distinguish the spiked peptides from the peptides originating from the sample. As preparation for nanoLC-MSMS analysis the peptides were purified by solid-phase extraction.

Results

Adjusting the “simple area based ratio estimation” using heavy isotopically labelled peptide standards, eliminated most of the previous observed deviation. By spiking the heavy isotopic labelled peptides in known concentrations, it should be possible to quantify the amounts of variants in the sample, corrected for differing MS response and possible sample matrix effect. However, recovery of the spiked internal standard is not the only challenge of an absolute quantification, which also depend on recovery of the many steps in the sample preparation procedure (solubilization/extraction, BCA protein quantification, digestion and SPE purification).

Conclusions

The precision, of a “simple chromatographic peak-area-ratio” describing the ratio between two protein variants in a sample, was remarkably improved by adjusting to internal standards in the form of heavy isotopic labelled peptides. Absolute quantification depends on many steps of the sample preparation procedure and each step require validation.

TP-69 / THINK POSITIVE: FINDING THE BEST MATRIX IN MALDI-TOF MS FOR CROCONAINE DYES INVESTIGATION

Cosima Damiana Calvano - Maria Capozzi - Angela Punzi - Gianluca Farinola - Francesco Palmisano - Tommaso Cataldi

University of Bari, Chemistry, Bari (Italy)

Keywords: croconaine, electron transfer, proton transfer, MALDI

Introduction

Croconaine dyes are appealing molecules based on the condensation of croconic acid and a reactive electron-donating aromatic or heterocyclic structures [1,2]. Fast investigation of these dyes was accomplished by matrix assisted laser desorption/ionization (MALDI) MS in positive ion mode after a systematic matrix selection. Since few or no information on the gas-phase ion energetics of croconaines are available, an empirical approach was followed.

Methods

All experiments were performed using a 5800 MALDI ToF/ToF analyzer in reflectron positive mode. Different matrices such as protonating 2,5-dihydroxybenzoic acid (DHB), and α -cyano-4-hydroxycinnamic acid (CHCA), electron-transfer (ET) secondary reaction matrices as trans-2-[3-(4-t-butyl-phenyl)-2-methyl-2-

propenylidene]malononitrile (DCTB) and 1,5-diaminonaphthalene (DAN) and basic as 9-aminoacridine (9AA) were examined.

Results

The investigated croconaines bring symmetrical substituents on the indolenine core including heteroatoms as fluorine and bromine or different aromatic structures. Typically, when dyes were analyzed by using proton transfer (PT) matrices as DHB and CHCA the base peak assigned to the protonated adduct $[M+H]^+$ and the concomitant presence of an odd electron molecular ion, $M^{+\bullet}$, was observed. This outcome suggests that gas-phase PT and electron transfer ET secondary reactions are competing processes when acidic matrices are employed. The presence of many signal contributions represents a drawback in terms of sensitivity and trouble when running tandem MS. Among the ET matrices, DAN provided the highest ionization yield leading to the specific formation of odd-electron molecular ion $M^{+\bullet}$ without other contributions. This allowed to select the precursor ion in a unique form for MS/MS analyses. Singular behavior was observed for bromine substituted croconaine when fragmented with DAN since bromine ion was detected as unusual Br^+ cation.

Conclusions

Protonating matrices generate a mix of odd molecular ions and protonated adducts alongside sodiated species. Among the ET matrices, DAN provided the highest ionization yield leading to the specific formation of odd-electron molecular ions $M^{+\bullet}$ with a negligible fragmentation. Finally, MALDI MS/MS of the radical charged molecular species provide useful structural information, thus making identification very straightforward for all croconaines investigated.

Novel Aspect

The analysis of croconaines by MALDI MS is presented for the first time. Tandem MS and DAN reveal common fragmentation patterns and the occurrence of uncommon positive ions.

References

1. Punzi A., Capozzi M. A. M., Fino V., et al., *J. Mater. Chem. C*, 4, 3138–3142 (2016).
2. Simard Y., Zebrowski Y., Haley D., *J. Org. Chem.* 65, 2236-38 (2000)

TP-70 / EFFICIENT PREPARATION OF ORGANOIMIDO DERIVATIVES OF LINDQVIST HEXAMOLYBDATE IN LEIDENFROST DROPLETS

Jie Cao - QianQian Wang

Key Laboratory of Cluster Science, Ministry of Education of China; Beijing Key Laboratory of Photoelectronic/Electrophotonic Conversion Materials; School of Chemistry and Chemical Engineering, Beijing Institute of Technology, Beijing, 100081, China

Keywords: Efficient Preparation; Organoimido functionalization; Lindqvist hexamolybdate; Leidenfrost droplets; Electrospray mass spectrometry

Introduction

Polyoxometalates (POMs) are discrete anionic metal-oxygen clusters which exhibit a great diversity of applications ranging from catalysis, medicine, etc. [1] POM-based inorganic/organic hybrids have drawn enormous attention in the past decades as exciting applications of these materials have been discovered in such diverse areas as light-emitting diodes, field-effect transistors, and solid-state lasers. Such hybrid materials may provide additional and/or enhanced functions and properties as a result of synergistic interactions between their inorganic and organic components.

Peng et al. reported a facile and efficient method to synthesize organoimido derivatives of POMs by using

dicyclohexylcarbodiimide (DCC) as a catalyst, which was implemented by refluxing the reaction mixture of aromatic amines and TBA₂[Mo₆O₁₉] in acetonitrile for 12 hours. [2] Leidenfrost droplets have been used as microreactors for the preparation of metal nanoparticles and for accelerated organic synthesis. [3] However, there is no report on the preparation of inorganic-organic hybrids which show promise in material science and homogeneous catalysis. Herein, we report preliminary result for the functionalization of POM clusters with aromatic amines, a model example of this type of reaction and its first implementation using the Leidenfrost effect to perform accelerated chemical synthesis in droplets. We show that interfacial effects are involved in acceleration and prepared macroscopic amounts of synthetic products within minutes. Compared to Peng's method, which was reported to be faster, more efficient, and the reaction conditions are much milder than earlier methods, Leidenfrost experiment is more competitive in which interfacial effects may play a pivotal role in the reaction acceleration.

Methods

Efficient Preparation of Organoimido Derivatives of Lindqvist hexamolybdate in Leidenfrost Droplets.

Results

Four reactions of Lindqvist hexamolybdate [Mo₆O₁₉]²⁻ with aromatic amines to generate organoimido derivatives were used as model reactions to explore the use of the Leidenfrost effect for accelerated chemical synthesis of organo-modified POMs in droplets. It is clear to see that observable amounts of products were formed over the time course of the experiment (240 s) for the Leidenfrost droplet experiment while in the parallel bulk-phase reactions no detectable product formation is evident.

Conclusions

We show that interfacial effects are involved in acceleration and prepared macroscopic amounts of synthetic products within minutes. Compared to Peng's method, which was reported to be faster, more efficient, and the reaction conditions are much milder than earlier methods, Leidenfrost experiment is more competitive in which interfacial effects may play a pivotal role in the reaction acceleration.

Novel Aspect

We report the first implementation using the Leidenfrost effect to perform accelerated synthesis of organo-functionalized POMs in droplets.

References

1. M. T. Pope, *Heteropoly and Isopoly Oxometalates*; Springer-Verlag: New York, 1983
2. *Polyoxometalates: from Platonic Solids to Anti-Retroviral Activity*; M. T. Pope, A. Müller, Eds. ; Kluwer Academic Publishers: Dordrecht, The Netherlands, 1994.
3. Y. G. Wei, B. B. Xu, C. L. Barnes, and Z. H. Peng, *J. Am. Chem. Soc.* 123, 4083-4084 (2001).
4. R. Abdelaziz, D. Disci-Zayed, M. K. Hedayati, J. H. Pöhls, A. U. Zillohu, B. Erkartal, V. S. K. Chakravadhanula, V. Duppel, L. Kienle, M. Elbahri, *Nat. Commun.* 4, 2400 (2013)
5. D. W. Lee, M. H. Jin, C. B. Lee, S. W. Lee, J. W. Park, D. Oh, J. C. Park, J. S. Park, *Chem. Mater.* 27, 5151-5160 (2015).

TP-71 / A NOVEL PTR-TOF REACTION CELL SUPERPOSING DC AND RF FIELDS LEADS TO TEN-FOLD IMPROVEMENTS IN SENSITIVITIES AND LODS

Luca Cappellin⁽¹⁾ - ***Felipe Lopez***⁽¹⁾ - ***Manuel Hutterli***⁽¹⁾ - ***Jordan E. Krechmer***⁽²⁾ - ***Sonja Klee***⁽¹⁾ - ***Benoit Plet***⁽¹⁾

TOFWERK AG, R&D, Thun (1) - Aerodyne Research Inc., R&D, Billerica (2)

Keywords: PTR-MS, time-of-flight, sensitivity, reaction cell, VOCUS

Introduction

The working principle of Proton Transfer Reaction - Mass Spectrometry (PTR-MS) is chemical ionization upon reactions between H_3O^+ ions and the analyte VOCs within a reaction cell. A DC field within the cell moves the primary and product ions towards the detector. In conventional reaction cell ion losses occur at the reaction chamber cell walls because of diffusion and scattering. This problem affects the instrumental sensitivities and limits of detection (LOD).

Methods

The solution involves adding oscillating RF fields on top of the linear field in the reaction cell. The ions are not only moved along the cell by the conventional DC field but are also simultaneously focused by the oscillating RF fields, so that more ions can reach the detector. The chemistry of ion-molecule reactions within the cell is similar to that of the old design, while the spatial distribution of primary and product ions is changed to achieve the extra focusing.

Results

The additional ion focusing which is reached in this way greatly improves the transmission of ions towards the detector region. The net gain in sensitivity exceeds one order of magnitude and previous boundaries in LODs are overcome thus providing unprecedented performances for a PTR-MS. The new reaction chamber has been coupled to state-of-the-art mass analyzers such as time of flight mass spectrometers (ToF-MS) reaching 15,000 ($m/\Delta m$ FWHM) in mass resolution. This new instrument has been denominated VOCUS PTR-TOF. In the case of xylene, the achieved sensitivity is 20000 cps/ppv and the LODs are <1 pptv and <10 pptv in 1 s and 1 min integration time, respectively. Another problem which is overcome by the new reaction cell is the dependence of calibration factors on the sample humidity. Calibration factors for VOCUS PTR-ToF are independent of sample humidity.

Conclusions

A new PTR-MS reaction cell design is introduced adding oscillating fields to the linear field. Ion losses due to scattering and diffusion are avoided, strongly improving the instrument performance. Limits of detection and sensitivity are improved by more than a factor of ten. The dependence of concentration calibration factors on sample humidity, which is a problem in conventional PTR-MS, is overcome.

Novel Aspect

The addition of an oscillating field to the PTR-MS reaction cell leads to a sensitivity gain of more than ten-fold.

TP-72 / DETERMINING ACETYLCHOLINESTERASE ACTIVITY USING PAPER SPRAY MASS SPECTROMETRY

Daniel Carmany⁽¹⁾ - ***Gabrielle Boyd***⁽²⁾ - ***Phillip Mach***⁽²⁾ - ***Elizabeth Dhummakupt***⁽³⁾ - ***Paul Demond***⁽²⁾ - ***Michael Busch***⁽²⁾ - ***Trevor Glaros***⁽⁴⁾

Excet Inc, ECBC, Gunpowder, MD (1) - Excet Inc, ECBC, APG, MD (2) - Nation Research Council, ECBC, APG, MD (3) - US Army, ECBC, APG, MD (4)

Keywords: Paper spray, acetylcholinesterase, chemical warfare agents, VX, enzyme

Introduction

Paper spray (PS) is an ambient ionization technique that has been used to detect and quantitate a wide range of small molecules from complex matrices. Recently, PS mass spectrometry was shown to be a rapid alternative for measuring the results of enzymatic reactions which took place in microcentrifuge tubes. Herein, we present an assay where the enzymatic reaction occurs on the PS substrate to monitor the activity of acetylcholinesterase (AChE) in blood to determine inhibition due to chemical warfare agent

exposure. This substrate assay is simple, rapid, and requires minimal sample volume.

Methods

The paper substrate required vapor application of silanes to generate a hydrophobic surface. The reaction was generated in a vial and spotted onto the surface of the PS cartridge. Due to the hydrophobicity of the substrate surface, the reaction solution formed a 'reaction sphere'. The loaded PS cartridges were then incubated at 37°C for 30 minutes. The cartridges were then analyzed on a Thermo Scientific Orbitrap Elite mass spectrometer for AChE activity.

Results

Enzymatic reactions were successfully performed on-substrate using PS-MS to determine AChE activity. AChE activity was determined using the synthetic substrate 1,1-dimethyl-4-acetylthiomethylpiperidinium (MATP), m/z of 202.0836, which is cleaved by AChE generating a resultant fragment at m/z of 160.0837. The change in ratio of these 2 ions demonstrates AChE activity, with corresponding increases in signal as additional substrate was reacted. It was determined that a 30 min incubation at 37°C resulted a log difference in the ratio. The nerve agent VX was used as a control for the reaction to identify AChE as the functioning enzyme, via direct inhibition. A method was then designed that allowed for not only the determination of AChE activity but to also identify the toxic inhibitor, VX, with tandem mass analysis. This method was used to analyze blood from animals exposed to VX. The results were compared with analysis of the same samples utilizing the standard Ellman's Assay. There was consensus between the methods.

Conclusions

PS-MS is suitable for determination of AChE activity in blood. The assay requires only MTAP as a reagent, and no secondary reaction is required. Using a hydrophobic PS cartridge, the reaction may be performed on-substrate with both high and low resolution instruments being utilized for analysis. The assay performed similarly to standard laboratory procedures when analyzing samples from exposed animals, even at a dose a log less than is required to observe symptoms.

Novel Aspect

An assay has been developed to measure AChE inhibition using PS-MS. The assay requires a single reagent and also allows for mass based identification of the toxicant.

Funding for this project is provided by the Defense Threat Reduction Agency – Joint Science and Technology Office for Chemical and Biological Defense

TP-73 / LC-HRMS BY MEANS OF UNTARGETED DATA-DEPENDENT ACQUISITION FOR THE CHARACTERIZATION OF 'SUPERFOODS'

Francesca Di Ottavio⁽¹⁾ - *Julia Gauglitz*⁽²⁾ - *Louis Felix Nothias*⁽²⁾ - *Morgan Panitchpakdi*⁽²⁾ - *Pieter Dorrestein*⁽²⁾ - *Dario Compagnone*⁽¹⁾ - *Manuel Sergi*⁽¹⁾

University, University of Teramo/ Faculty of Bioscience, Teramo (1) - UCSD, University of California San Diego/ Skaggs School of Pharmacy, San Diego (2)

Keywords: Mass-Spectrometry, Food Analysis, super-food, GNPS, in-silico annotation

Introduction

We provide a workflow to identify characteristic signatures of superfoods, a poorly defined marketing term referring to food with high levels of nutrients that may benefit health.

A set of 500 plant-based and superfood samples were analyzed by using untargeted UPLC-QTOF-MS/MS and new computational tool. We identified distinct molecular signatures including the diversity of

polyphenols, flavonoids and microbial peptides that can be used to categorize superfoods.

Methods

Samples were extracted and chromatographically separated by UPLC-QTOF system and C18 column separation. Spectra were acquired in positive data dependent acquisition and the five most abundant ions were selected for fragmentation spectra acquisition. Greater than 500 plant-based samples were analyzed by UPLC-MS/MS in DDA and processed on the GNPS web platform [1].

Results

Untargeted MS analysis was used to obtain an overview of all features, focusing on compounds with anti-inflammatory, antioxidant or antimicrobial properties across a broad range of plant-based foods. Spectral library search was performed against public and commercial libraries (GNPS, MassBank, Metlin, NIST), and molecular networks were generated to expand spectral annotations. In addition, computational mass spectrometry strategies were employed to annotate bioactive compounds such as polyphenols, flavonoids and phenolic acids. Those include: (1) the automatic network annotation propagation of in silico annotations, (2) the dereplication of bioactive peptides using the Dereplicator workflow [2], (3) the annotation of sub-structures using MS2LDA [3], (4) the computation of the molecular formula using Sirius and CSI:Finger [4] and (5) the prediction of the chemical class using Canopus fingerprint.

Conclusions

Unique microbial peptides with antimicrobial, anti-inflammatory or antioxidant properties and putative unknown cyclolinopeptide variants were found in some superfoods (score>15, p values< 1E-11, FDR< 1%). The presence of unique compounds, as well as the proportions of polyphenols, flavonoids and bioactive peptides between superfoods and common fruit and vegetables is proposed to define criteria underlying the chemical characteristics of superfoods.

Novel Aspect

We are presenting criteria to classify superfoods based on the number and the diversity of bioactive compounds.

References

1. Wang M. et al., Nature Biotechnology, 34, 828-837 (2016).
2. Mohimani H., Gurevich A., Mikheenko A., Garg N., Nothias L. F., Ninomiya A., Takada K., Dorrestein P. C., Pevzner P. A., Nature Chemical Biology, 13, 30-37 (2017).
3. Van Der Hooft J. J., Wandy J., Barrett M. P., Burgess K. E., Rogers S., Proceedings of the National Academy of Science of the United State of America, 113 (48), 13738- 13743 (2016).
4. Duhrkop K., Scheubert K., Bocker S., Metabolites, 3 (2), 506-516 (2013)

TP-74 / HIGH-RESOLUTION MASS SPECTROMETRY IN COMBINATION WITH IN SILICO PREDICTION MODELS TO IMPROVE ACCURACY FOR COMPOUND IDENTIFICATION

Eric Dossin - Pierrick Diana - Anneke Glabasnia - Mark Bentley - Philippe Alexandre Guy

Philip Morris Product S. A., Product research, Neuchâtel

Keywords: compound identification, high-resolution mass spectrometry, complex mixture analysis, linear retention index prediction

Introduction

Chemical characterization remains a major hurdle for analytical chemists when dealing with non-targeted

screening analyses. This work focuses upon the identification of unknown volatile and semi-volatile compounds using gas chromatography coupled to high-resolution mass spectrometry (GC-HRMS) operating in electron ionization (EI) acquisition mode.

Methods

552 reference compounds, including n-alkanes used as chemical markers, were analyzed by GC-HRMS. Accurate mass spectra (MS) were registered in a Personal Compound Database Library, which enabled targeted screening for these compounds. These compounds were randomly split into training (n=401) and test (n=151) sets, and linear retention index (LRI) values were used to build LRI prediction models [1].

Results

These models were used to predict LRI values for several thousands of chemicals reported to be present in tobacco and tobacco aerosol, plus a range of specific flavor compounds[2]. These predicted LRI values were used in conjunction with EI nominal MS available from commercial libraries to enhance the confidence in structural elucidation.

These models were also used to predict LRI values for hits proposed from the comparison of deconvoluted spectra with commercial MS libraries. The similarity between calculated and predicted LRI values added an additional layer of confidence for compound identification. A final score was calculated as a combination of MS similarity score and the accuracy for experimental versus predicted LRI values.

Recently, a set of 350 reference compounds were analyzed and used to validate the accuracy of the models and highlight any limitations.

Conclusions

The combination of accurate LRI prediction and MS similarity improved the ranking of compound proposals, thus demonstrating that computational LRI prediction facilitates the evaluation of non-targeted data. This workflow has been used to improve the confidence in structural elucidation for the analysis of complex matrices.

Novel Aspect

Combination of two LRI prediction models with EI acquisition mode for identification of unknown chemicals.

References

1. Dossin E., Martin E., Diana P., Castellon A., Monge A., Pospisil P., Bentley M., Guy P. A. *Anal. Chem.* 88,15, 7539-7547 (2016).
2. Martin E., Monge A., Duret J., Gualandi F., Peitsch M. C., Pospisil P. J. *Cheminform.* 4, 1-11 (2012).

TP-75 / HIGH-THROUGHPUT QUAN/QUAL LC-MS ANALYSIS: WHERE IS THE SWEET SPOT?

Anne-Charlotte Dubbelman⁽¹⁾ - ***Filip Cuyckens***⁽²⁾ - ***Lieve Dillen***⁽²⁾ - ***Rob J. Vreeken***⁽²⁾ - ***Thomas Hankemeier***⁽¹⁾

Leiden University, Leiden Academic Centre for Drug Research/Systems Biomedicine and Pharmacology, Leiden (1) - Janssen R&D, Pharmacokinetics, Dynamics and Metabolism, Beerse (2)

Keywords: Quan/Qual; UHPLC-QTOF-MS; high-throughput; metabolomics; drug metabolism

Introduction

High-throughput simultaneous quantitative and qualitative (Quan/Qual) analysis is attractive to combine targeted with non-targeted analysis, *e.g.* in pharmacometabolomics and drug metabolism studies. But, when using a Q-ToF, where is the 'sweet spot', viz. an optimal balance between quantitative and qualitative data and optimal use of instrumental capabilities? That is what we investigated in this study.

Methods

An ultra-high performance liquid chromatography (UHPLC) method with high chromatographic resolution (peak capacity of 266 with a gradient time of 2.25 min) and a Synapt G2S were used to investigate and optimize mass spectrometric parameters for Quan/Qual analysis. The optimized method was applied to a rat metabolomics study investigating the effect of the fasting state and the administration of a dosing vehicle on the rat plasma metabolic profile.

Results

The scan mode and scan rate proved to be highly important parameters for high-throughput Quan/Qual analysis. The amount of qualitative information that a method can provide correlated negatively with its quantitative performance (accuracy, precision, sensitivity, linear dynamic range). The MSE scan mode with a short scan time of 30 ms provided the optimal balance. The resulting 4.25-min Quan/Qual analysis method enabled quantification with accuracy and precision values $\leq 20\%$ at the lowest quality control (QC) level and $\leq 15\%$ at higher QC levels for 16 out of 19 tested analytes. It provided both parent m/z values and fragmentation spectra for compound identification with limited loss of chromatographic resolution. It revealed biologically relevant metabolites in its application to the rat metabolomics study.

Conclusions

High-throughput Quan/Qual method development requires balancing between the amount of qualitative data, the quality of the quantitative data and the analysis time. We will share recommendations for MS resolution, scan mode, scan rate and smoothing in Quan/Qual method development and analysis.

Novel Aspect

This work highlights the importance of optimizing high-resolution mass spectrometric parameters for Quan/Qual analysis, opposing the general idea of their insignificance.

TP-76 / UHPLC-HRMS METHOD FOR TARGETED ANALYSIS OF CYANOTOXINS IN FRESH WATER

Daria Filatova⁽¹⁾ - *Massimo Picardo*⁽¹⁾ - *Oscar Núñez*⁽²⁾ - *Marinella Farré*⁽³⁾

Institute of Environmental Assessment and Water Research, Spanish National Research Council, Department of chemical engineering and analytical chemistry, University of Barcelona, Barcelona (1) - Serra Hunter Fellow, Generalitat de Catalunya, Department of chemical engineering and analytical chemistry, University of Barcelona, Barcelona (2) - Institute of Environmental Assessment and Water Research, Spanish National Research Council, Barcelona (3)

Keywords: cyanotoxins, microcystins, high-resolution mass spectrometry

Introduction

Cyanobacteria are the components of regular microbial succession in periphyton formation, however, some species are known to produce toxic secondary metabolites (cyanotoxins), which vary in structure and harmful properties. During the past decades, an increase of cyanobacterial blooms was noticed in freshwater bodies. To monitor cyanotoxins' levels and prevent both human poisoning and wildlife damage, suitable analytical methods need to be developed.

Methods

For sample pretreatment, SPE using 2 cartridges in sequence for multiple toxins has been optimized and employed. The chromatographic separation was achieved using a C18 analytical column (150x2.1 mm, 2 μ m) with acidified acetonitrile and water as mobile phase. The chromatographic separation was coupled to a Q-Exactive Orbitrap instrument (Thermo Fisher Scientific) and the MS data was acquired in both full scan and parallel-reaction monitoring modes (PRM).

Results

Ultra-high performance liquid chromatography coupled to high-resolution mass spectrometry (UHPLC-HRMS) was developed. This method provides targeted analysis of cyanotoxins (microcystins, nodularin, cylindrospermopsin, and anatoxin-a) with instrumental limits of detection between 0.1 ng/L (4 toxins) and 1 µg/L (2 toxins). Application of PRM gave an opportunity to identify the most abundant fragments for each toxin. Additionally, HRMS provides an assessment to the potential presence of transformation products and other non-targeted toxins in the samples. The developed method was applied for the study and characterization of cyanotoxins concentrations in Catalonia freshwater reservoirs.

Conclusions

The developed method can be applied for identification and quantification of cyanotoxins of different chemical classes.

Novel Aspect

A sensitive, fast and robust method for the analysis of cyanotoxins in freshwaters based on UHPLC-HRMS was developed.

TP-77 / DIRECT GAS ANALYSIS USING A MULTI-TURN HIGH-RESOLUTION MASS SPECTROMETER

Masahiro Hashimoto⁽¹⁾ – **Marcel J. van der Schans**⁽²⁾ - **Jeroen A. Van der Meer**⁽²⁾ - **Bram van der Meer**⁽³⁾

JEOL (EUROPE) SAS, European Application Group, Croissy-sur-Seine (1) - TNO, Defence, Safety and Security, Rijswijk (2) – JEOL (EUROPE) B. V., Sales Executive Science and Measurement Instruments, Nieuw-Vennep, (3)

Keywords: Gas analysis, High resolution, High mass accuracy, Multi-turn

Introduction

Small molecule gases such as CO, NO_x and SO₂ are analyzed by using a portable gas analyzer or a GC-MS. However, portable gas analyzers are only able to measure specific gases, and GC-MS is able to measure a variety of gases but requires special columns to separate small molecule gases. In this study, we measured complex gas mixtures by using the multi-turn HR-TOFMS without special columns.

Methods

We measured gas mixtures by using the multi-turn HR-TOFMS with three different sampling methods - gas mixing system, headspace and a tedlar bag. For the gas mixing system, a sample was mixed that contained benzene (1250 ppm) and arsine (40 ppm). The gases were introduced by using the MS vacuum system to draw the sample directly into the source with a short deactivated capillary column.

Results

Target gases were observed, and their compounds were separated by the multi-turn HR-TOFMS perfectly. The mass error for the peaks were within 3.1 mDa. In the case of gas mixing system, the *m/z* difference of benzene and arsine is 0.1 Da (benzene, C₆H₆, M⁺· 78.04640; arsine, AsH₃, M⁺· 77.94452). The results showed that the overloaded peak for benzene (M⁺·) and the smaller peak for arsine (M⁺·) were perfectly separated and detected by the multi-turn HR-TOF MS. In addition, the arsine (M⁺·) mass error was only 0.5 mDa.

Conclusions

The multi-turn HR-TOFMS can separate complex gas mixtures by using high mass resolution and can identify elemental compositions by using the system's high mass accuracy. This means that the multi-turn HR-TOFMS is a powerful tool for analyzing complex gas mixtures.

Novel Aspect

Complex gas mixtures analysis by using a multi-turn HR-TOFMS.

TP-78 / DEVELOPMENT OF AN LC-HRMS METHOD TO AID IN IDENTIFICATION AND CONFIRMATION OF CHLORINATED AND BROMINATED CONTAMINANTS

Alin Ionas⁽¹⁾ - *Ana Miralles-Marco*⁽¹⁾ - *Garry Codling*⁽¹⁾ - *Simon Vojta*⁽²⁾ - *Lisa Melymuk*⁽¹⁾ - *Jana Klánová*⁽¹⁾

Research Centre for Toxic Compounds in the Environment (RECETOX), Masaryk University, Brno (1) - Graduate School of Oceanography, University of Rhode Island, Kingston (2)

Keywords: halogenated chemicals, LC-HRMS, compound identification, Orbitrap Fusion, Data-dependent analysis

Introduction

In human and environmental exposure, many of the compounds of concern for their environmental persistence and health risks are halogenated. The twelve initial persistent organic pollutants (POPs) included in the Stockholm Convention are all chlorinated. Many of the new POPs added subsequently are brominated, such as hexabromobiphenyl, (HBB), hexabromocyclododecane (HBCD), polybrominated diphenyl ethers (PBDEs, from the Tetrabrominated to Heptabrominated), or fluorinated, such as perfluorooctane sulfonic acid (PFOS), its salts and perfluorooctane sulfonyl fluoride (POSF) [1]. Furthermore, most chemicals listed in studies evaluating the potential of compounds to meet the criteria of being a POP, (persistence, bioaccumulative and ability to undergo long-range transport) are also halogenated [2,3,4].

For this reason, in the past few years there have been several studies trying to identify and highlight previously unknown halogenated chemicals of both natural and anthropogenic sources such as in stream sediments from southern Ontario [5], Atlantic dolphin blubber [6], and *Prymnesiumparvum* haptophytes [7] and others.

In general, for this type of studies, the samples are usually acquired in full scan or data-independent analysis mode. This is followed by a lengthy data processing procedure, typically involving three steps:

- 1) Target screening / quantitative analysis step, where information about the levels of analytes present are confirmed with standards
- 2) Suspect screening step, where compounds of interest (for which standards are not available) are monitored
- 3) Non-target screening step, where the identification of other compounds is attempted.

For the confirmation of chemicals tentatively identified, tandem mass spectrometry information is necessary. To obtain it, an additional targeted MS_n analysis is often required, even if the initial acquisition mode was data-independent analysis, if the intensity of the peak of interest is not high enough to be selected for fragmentation.

Consequently, the aim of this study is to develop an instrumental method designed to speed up and facilitate later data analysis by:

- 1) Incorporating an element of target and suspect screening in the acquisition method itself
- 2) Highlighting chemicals containing more than three atoms of chlorine and bromine, which are more likely to have an adverse effect on human health
- 3) Obtain MS₂ scans of the halogenated compounds of interest to use in their identification, without requiring subsequent instrument time.

Methods

The instrument employed in the present study was a Shimadzu Nexera X2 UHPLC coupled to an Orbitrap Fusion MS detector (Thermo Fisher Scientific, Waltham, MA, USA). The detector was run at a resolution of 120,000 and was equipped with an atmospheric pressure chemical ionization (APCI) source. As source parameters, the negative ion discharge current was set at 10 μ A and the positive at 4 μ A. The sheath gas was set at 25 arbitrary units (equivalent to 3.17 L/min), and the auxiliary gas at 5 arbitrary units (corresponding to 5.08 L/min). The sweep gas was turned off. The ion transfer tube temperature was set at 250°C and the vaporizer temperature at 325°C.

A volume of 5 μL of extract was injected, and separation was achieved using a Waters Acquity BEH C18 column (100 mm \times 2.1 mm i.d., 1.7 μm particle size), fitted with a 5 mm matching Vanguard pre-column, with the same characteristics. The mobile phases employed were A: milliQ water and B: methanol, both with 5 mM ammonium formate added. The flow rate was 0.3 mL/min and a linear gradient from 5% to 100% methanol in 8 min, followed by a 2 min hold before returning to the initial conditions for 4 min. Standards were injected both in negative and positive polarity. For testing and development of the method, standards of brominated and chlorinated chemicals – mainly flame retardants – were employed. The “BDE-CM” mixture of PBDE congeners (28, 47, 99, 100, 153, 183 and 209) was purchased from AccuStandard (New Haven, CT, USA), Hexabromocyclododecane (HBCD) from Wellington Laboratories (Guelph, ON, Canada) and three mixtures of novel FRs (From the INTERFLAB inter laboratory study [8]) were kindly donated by Wellington Laboratories. The mixtures contained 26 FRs in total: 2,4,6-Tribromophenyl allyl ether (ATE), 2,3,5,6-Tetrabromo-p-xylene (p-TBX), 2-bromoallyl-2,4,6-tribromophenyl ether (BATE), 1,2,3,4,5-pentabromobenzene (PBBz), Tetrabromo-o-chlorotoluene (TBCT), Pentabromotoluene (PBT), Pentabromoethylbenzene (PBEB), Tris(1,3-dichloroisopropyl) phosphate (TDCPP), 2,4,6-Tribromophenyl 2,3-dibromopropyl ether (DPTE), Hexabromobenzene (HBBz), Tris(2-butoxyethyl) phosphate (TBEP), Pentabromobenzyl acrylate (PBBA), 2-Ethylhexyl 2,3,4,5-tetrabromobenzoate (EHTBB), Tetrabromobisphenol A (TBBPA), 1,2-Bis(2,4,6-tribromophenoxy)ethane (BTBPE), Bis(2-ethylhexyl)tetrabromophthalate (BEHTBP), syn-Dechlorane Plus (s-DP), alpha-Dechlorane Plus (a-DP), Octabromotrimethylphenylindane (OBIND), Decabromodiphenyl ethane (DBDPE), α -Tetrabromoethylcyclohexane (aTBECH), β -Tetrabromoethylcyclohexane (bTBECH), α -1,2,5,6-Tetrabromocyclooctane (aTBCO), β -1,2,5,6-Tetrabromocyclooctane (bTBCO), 5,6-Dibromo-1,10,11,12,13,13-hexachloro-11-tricyclo[8.2.1.0^{2,9}]tridecene (HCDBCO) and Tris(2,3-dibromopropyl)isocyanurate (T23BPIC; abbreviations are those used by Wellington Labs) considered as current use compounds. All acronyms are as described in the certificates of analysis from the manufacturers and reference [8].

Results

To achieve the aims described above, the detector needed to be fast enough to cycle between full MS acquisition and data-dependent MS₂ acquisition and the cycle time low enough to permit at least two cycles of MS₂ acquisition per peak. A cycle time of 0.8 seconds was chosen for this purpose.

Several filters were set after the initial Full MS scan, serving to select only the desired chlorinated and brominated ions to the latter type of acquisition.

The most relevant filters used to select just the ions with >3 atoms of Cl and/or Br for fragmentation and acquisition of MS₂ spectra were: 1) “Targeted Mass Difference” and 2) “Targeted Isotopic Ratio”. In the first filter, it was specified to select pairs of ions with a mass difference of 3.9941 Da 3.9959 for chlorinated and bromine-containing substances respectively. These values were chosen instead of the typical A/A+2 pair because of having the benefit of decreasing the number of false positives of non-chlorinated or brominated ions being selected for MS₂.

For the second filter, “Targeted Isotopic Ratio” the input was the typical ratios of isotopomers in a halogenated cluster of 3-14 atoms of Cl or Br. The ratios were calculated theoretically as described in a previous publication [9]. Some tests have also been done using an average value for these ratios and higher tolerance values, but this setup was less than ideal, as it generated higher numbers of false positives.

An exclusion list was defined by monitoring the most abundant ions coming from the mobile phase and the LC system. The list was uploaded into the “Targeted Mass Exclusion” filter, to prevent these ions from being selected for MS₂ fragmentation.

An “Intensity” filter was used as well, to eliminate ions close to background levels (filter set at 5E3).

The “Precursor Selection Range” was set at 250-1,200, to eliminate the 79/81 cluster and the 159.83 characteristic to many brominated analytes and thus highly unselective.

The “Dynamic Exclusion” filter had the purpose to ensure that ions from more than one halogenated cluster was selected for the data-dependent MS/MS scan, if more than one such cluster is present at the same retention time. This is accomplished by preventing an ion from triggering a second data-dependent scan, for 20 seconds after it had triggered the first one.

To increase the likelihood of obtaining high quality MS/MS data, the “Apex Detection” filter was used, to

trigger the subsequent data-dependent MS/MS scans after the apex of the peak, when the intensity of the clusters of interest is at its highest. For this purpose, an expected peak FWHM of 3 seconds was found to be adequate, with a desired apex window of 15%.

A second data-dependent MS/MS was added, preceded by a "Targeted Mass Trigger". This particular scan event was triggered whenever one of the masses on the list was detected. The ions observed in APCI negative mode for the analytes specified in the Materials and methods section were added to this list (targeted screening) and also accurate masses for other relevant analytes from the literature [9,10], such as 2,4,6-tris(2,4,6-tribromophenoxy)-1,3,5-triazine (TTBP-TAZ), for which standards were not available (suspect screening).

At the present stage of the method development process, the method correctly selects chlorinated and brominated clusters for data-dependent MS² in over half of the cases examined. This can be easily improved by increasing the tolerances of parameters such as the isotope ratio, targeted mass difference, decreasing dynamic exclusion time and adjusting the parameters apex detection (or removing it entirely). However, doing so increases the number of false positives obtained. For example, for BEHTBP, the full scan MS spectrum was obtained, the most abundant isotopomer from the main cluster was automatically selected for the first data-dependent MS² scan event and the same one was successfully selected as being on the list of the targeted mass trigger. But the same cannot be said for the second, less abundant isotopic cluster from the mass spectrum.

Conclusions

As a conclusion, the method is accomplishing most of what it was intended for, but more fine tuning is needed for it to be fully suited for its intended purpose. Further studies on highly complex samples with this method are planned, after an additional tuning using indoor dust and waste water stream samples.

Novel Aspect

To the authors' knowledge, this is the first study describing a method which provides targeted and suspect screening information, highlights halogenated chemicals in the chromatograms and provides MS² information for the identification and confirmation of halogenated analytes, all in one run.

References

1. <http://chm.pops.int/TheConvention/ThePOPs>
2. Scheringer M, Stempel S, Hukari Sn, Ng CA, Blepp M, Hungerbuhler K (2012) Atmospheric Pollution Research, 3, 383–391.
3. Brown TN, Wania F (2008) Environmental Science & Technology, 42, 5202–5209.
4. Howard PH, Muir DC (2010) Environmental Science & Technology, 44, 2277–2285.
5. Pena-Abaurrea M, Jobst KJ, Shen L, McCrindle R, Helm PA, Reiner EJ (2014) Environmental Science & Technology, 48, 9591–9599.
6. Hoh E, Dodder NG, Lehotay SJ, Pangallo KC, Reddy CM, Maruya K (2012) Environmental Science & Technology, 46, 8001–8.
7. Andersen AJC, Hansen PJ, Jørgensen K, Nielsen KF (2016) Analytical Chemistry, 88, 12461–12469.
8. Melymuk L, Goosey E, Riddell N, Diamond ML (2015) Analytical and Bioanalytical Chemistry, 407, 6759–6769.
9. Ionas AC, Ballesteros-Gómez A, Leonards P. E. G, Covaci A (2015) Journal of Mass Spectrometry, 50, 1031–1038.
10. Ballesteros-Gómez A, de Boer J, Leonards PEG (2013) Analytical Chemistry, 85, 9572–9580.

TP-79 / COMBINATION OF DIFFERENT IONIZATION SOURCES EXPANDS THE ANALYSIS OF COMPLEX MIXTURES

Aikaterini Kondyli - Wolfgang Schrader

Max Planck Institut für Kohlenforschung, Mass Spectrometry, Mülheim an der Ruhr

Keywords: atmospheric pressure ionization, crude oil, high resolution mass spectrometry

Introduction

Despite the great variety of ionization sources available, there is no single analytical technique for fully analyze a crude oil. Ion discrimination and ion suppression are limiting the compound detection [1,2]. While ESI is the method of choice for the detection of polar compounds, APPI and APCI favors the non-polar and low-to-medium polar compounds respectively. Until now, there is few information on the combination of different sources.

Methods

A light crude oil was dissolved and diluted in toluene/methanol (1:1) at a final concentration of 250 ppm. Mass spectrometric analysis was performed on an Orbitrap Elite. The spectra were collected in SIM mode with a mass resolution of 480000 (at m/z 400, 1.5s transient). ESI positive, APPI and APCI as well as the combination of ESI with APPI, or APPI with APCI were selected as ionization methods. For APPI, a Krypton VUV lamp emitting radiation of 10 and 10.6 eV was used.

Results

The results obtained with different ionization methods show clearly that under ESI the ionization of nitrogen containing compounds is favored, while under APPI hydrocarbon species as well as oxygen containing compounds in both protonated and radical form were efficiently ionized. The combination of these two methods acted as an intermediate, allowing the ionization of non-polar to polar compounds. An overall increase was observed mainly in the O1 and O2 classes compared to the results obtained only from electrospray, while the classes of OS1 and S1 showed the highest population abundance after the combination of these methods. On the other hand, in case of APCI, the hydrocarbon species as well as the N1, O1, O2, OS1 and S1 classes showed a higher abundance in their protonated form compared to the APPI. After their combination, an increase of the radical species in all the classes was observed.

Conclusions

Although each ionization method favors the detection of specific types of compounds, the combination of these methods gives the opportunity to utilize all the benefits of the individual methods. The use of different ionization methods provides complementary information about the crude oil constituents. To obtain information on a molecular level, high resolution mass spectrometry is the method of choice for the analysis of such complex mixtures.

Novel Aspect

The effect of combining different ionization methods together for the analysis of complex mixtures is investigated on a molecular level by HRMS.

References

1. A. Gaspar, E. Zellermann, S. Lababidi, J. Reece, W. Schrader. Anal. Chem. 84, 5257-5267, (2012).
2. Cho, Y., et al. Energy & Fuels 26(5), 2558-2565,(2012).

TP-80 / ABERRATION-BASED NUMERICAL OPTIMIZATION OF ION MIRRORS

Igor Kopaev⁽¹⁾ - **Dmitry Grinfeld**⁽²⁾ - **Mikhail Monastyrskiy**⁽³⁾ – **Alexander Makarov**⁽²⁾

General physics Institute, T8, Krasnobogatyrskaya street, 44 Moscow, Russian Federation (1) - ThermoFisher Scientific, Global Research LSMS, Bremen, Germany (2) - General Physics Institute, Optoelectronics Departement, Moscow, Russian Federation (3)

Keywords: Mathematical modeling, multireflectron geometry optimization

Introduction

Synthesis and optimization of a multireflectron system based on ion mirrors

Mass-spectrometry systems based on ion mirrors have been introduced in 1973 [1] and now represent a perspective class of high-resolution MS equipment [2, 3]. The aim of this work is computer-aided optimization of a planar-symmetric systems of meshless ion mirrors with the use of different numerical approaches.

Methods

Numerical approaches to optimization of planar-symmetric meshless ion mirrors

A part of the Matlab PDE Tollbox program code oriented to solve the Laplace equation was updated and tested. One of the main features of the newly developed code is a possibility of performing numerical optimization as applied to ion mirror simulation based on aberration theory.

Results

The classes of ion mirror geometries satisfying the prescribed requirements have been found.

The possibility of further minimization of the aberrations of interest has been shown.

The alternative geometries have been studied and simplified.

Conclusions

As a result of optimization calculations, a new geometry of the "stairs" type has been designed, which represents a compromise between different geometries based on rectangular electrodes.

On the other hand, it is shown that the geometries with rectangular electrodes possess fairly small aberrations of interest and are the simplest among the geometries explored

Novel Aspect

Computational optimization algorithms based on aberration theory to synthesize ion multireflectron systems with prescribed characteristics have been implemented.

References

1. Mamyrin, B., et al., Zh. Eksp. Teor. Fiz., 64, 82 (1973).
2. Grinfeld, D., et al., European Journal of Mass Spectrometry, 20(2), 131 (2014).
3. Skoblin, M., et al., Quantum Electronics, 45(12), 1171 (2015).

TP-81 / APPLICATION OF ETHCD FOR LEUCINE/ISOLEUCINE DISCRIMINATION IN PEPTIDES

Sergey Kovalev - Miriam Tolpina - Sergey Zhokhov - Tatiana Samguina - Albert Lebedev

Moscow State University, Chemistry Department, Moscow

Keywords: ETHcD, Orbitrap, peptidomics, leucine/isoleucine discrimination

Introduction

For many years Edman degradation was a golden standard for peptide sequencing. However since 1990s mass spectrometry (MS) has become an effective alternative with significant advantages like possibility to uncover non-standard amino acids and post translational modification. One problem for MS in this field is isomeric leucine / isoleucine discrimination. The aim of our work is to develop reliable and efficient method exploiting ETD to overcome this difficulty.

Methods

Experiments were carried out on Orbitrap Elite mass spectrometer with electrospray ionization source. Tryptic peptides of the human serum albumin, bacterial protein gp188, natural peptides isolated from skin secretion of *Rana ridibunda* frog and several synthetic peptides were introduced into the ion source using a

syringe pump. Mass spectra were measured in a positive ionization mode in Orbitrap mass analyzer with 120,000 resolving power.

Results

As a starting point we improved MS3 sequential ETD-HCD activation developed by our group earlier. To achieve maximum spectral quality and fragmentation efficiency all parameters of ion collection and fragmentation energy in MS3 ETD-HCD mode with maximum possible ion isolation width (600 Da) for peptides were optimized. The resulting ETD-HCD method, with "broad-band" activation made it possible to identify structure of 83 out of 93 leucine / isoleucine residues. Interfering process of radical site migration may be suppressed by a careful choice of the lowest possible energy for w-ion formation. The main disadvantage of "broad-band" activation is low efficiency of z-type ion formation for long peptides (with m/z more than 2000 Da). In the advanced version ETHcD was used to accelerate data acquisition. Besides it improved the spectra quality and had a greater sensitivity. Application of new method to a non-separated mixture of skin secretion peptides of frog *Rana ridibunda* allowed correctly identifying 60 of 66 Leu / Ile residues.

Conclusions

New ETHcD method exhibiting nearly 90% efficiency was developed for reliable Leu/Ile discrimination in peptides. Aspects of ion optics tuning were studied and ETD, ETcaD, and ETHcD modes compared. Radical site migration interfering with Leu/Ile discrimination can be suppressed by careful choice of minimal NCE. The present version of the developed approach may constitute a basis for the creation of automatic sequencing protocol.

Novel Aspect

The developed method is advantageous compared to existing ones. It has a potential for automation and coupling to LC to become a routine proteomics protocol.

TP-82 / A NEW LC-MS/MS METHOD FOR DETERMINATION OF APIXABAN IN TABLETS

Murat Kozanli - Nafiz Öncü Can

Anadolu University, Faculty of Pharmacy, Department of Analytical Chemistry, Eskisehir

Keywords: apixaban, monolithic silica column, LC-ESI-MS/MS, pharmaceutical application

Introduction

Apixaban (APX), a selective, reversible and direct inhibitor of factor Xa, is in use to reduce the risk of stroke and systemic embolism in patients with non-valvular atrial fibrillation [1]. In spite of some published methods for biological samples [2,3], there are very few methods on its determination in pharmaceutical preparations. The aim of this study was to develop a validated LC-MS/MS method for determination of APX in tablets.

Methods

Analyses were performed using a Shimadzu LCMS-8040 (Japan) series LC-ESI-MS/MS instrument. APX and irbesartan (internal standard; IS) were separated on a monolithic high-resolution silica column using isocratic elution. The detection was performed in positive ion mode and the precursor-to-product ion transitions were monitored in MRM mode at 460.15→443.20 m/z for APX and 429.15→207.15 m/z for IS. The developed method was validated as per ICH Q2(R1) guideline.

Results

For APX, linearity could be demonstrated over the concentration range of 2–20 $\mu\text{g/mL}$ ($R^2 > 0.996$), and the

calculated limit of quantification were about 1 µg/mL. results for inter- and intra-day assay precision and trueness were obtained using internal quality control samples and remained within the acceptance criterion. External quality control samples were measured at the specified nominal values with inter- and intra-day precisions. Matrix effects were fully compensated by co-eluting IS, which in turn did not relevantly influence ionization efficiency.

Conclusions

The method development was focused on the optimization of both LC and MS/MS conditions to obtain the best possible sensitivity and future applicability on regular LC-PDA systems. The proposed method was successfully applied on the real samples, possessing very high chromatographic efficiency.

Novel Aspect

The utilization of a monolithic high-resolution silica column, lower quantitation limits and use of IS are novel aspects of the study.

References

1. Landge, S. B., Jadhav, S. A., Dahale, S. B., Solanki, P. V., Bembalkar, S. R., Mathad, V. T., American Journal of Analytical Chemistry, 6, 539-550 (2015)
2. Schmitz, E. M., Boonen, K., van den Heuvel, D. J., van Dongen, J. L., Schellings, M. W., Emmen, J. M., van der Graaf, F., Brunsveld, L., van de Kerkhof, D., Journal of Thrombosis and Haemostasis: JTH, 12, 1636-1646 (2014).
3. Gous, T., Couchman, L., Patel, J. P., Paradzai, C., Arya, R., Flanagan, R. J., Therapeutic Drug Monitoring, 36, 597-605(2014).

TP-83 / MULTIPLEXED OPERATION OF HIGH RESOLUTION ORTHOGONAL MULTI-REFLECTING TIME-OF-FLIGHT MASS SPECTROMETER

Boris Kozlov⁽¹⁾ - ***Jeffrey Brown***⁽¹⁾ - ***Viatcheslav Artaev***⁽²⁾

Waters Corporation, Research Department, Wilmslow, United Kingdom (1) - LECO, St Joseph (2)

Keywords: multi-reflecting TOF MS, duty cycle, encoded frequent pulsing, folded flight path, artifact-free inverse transform

Introduction

Multi-reflecting TOF MS incorporate extended ion flight paths to increase resolution, however this is at the expense of duty cycle. Upstream trapping to accumulate ions results in loss of mass accuracy. Alternatively, multiplexing of TOF spectra has been proposed by numerous researchers [1,2,3]. Here we describe the implementation of Encoded Frequent Pushing (EFP) [4, 5], where a repeated sequence of unique time interval pusher pulses is applied.

Methods

In contrast to Hadamard transform methods, EFP involves a pattern of TOF start pulses where it's shifted replica does not coincide at more than one point (with a precision comparable to TOF peak width). Decoding by robust statistical verification of all data points produces no artefacts. The main parameter of the decoding algorithm is the ion population of multiplexed spectrum.

Data were acquired using a MR-TOF MS experimental system having a 46m long Folded Flight Path.

Results

The 2 ms long time-of-flight spectra were acquired by a 1.8 Giga-samples/sec 12-bit ADC (X6 GSPS Innovative Integration), with on the fly decoding performed using computational resources of the GPU card (NVIDIA). They demonstrated two orders improvement of duty cycle and hence sensitivity with no loss of resolution and with dynamic range reaching 5 orders of magnitude in a single ESI-spectrum.

There was no degradation of the mass accuracy when using EFP recorded over the mass to charge range 300 to 3000 (using 2-points for internal calibration). In fact, for weaker peaks within the data, inclusion of more ions helped to improve the accuracy of the mass measurement.

We also found that EFP does not introduce deviations to the relative intensities of various spectral peaks, except when the signal was at the minimal detectable level.

Conclusions

The method is artefact free and provides duty cycles comparable with conventional, non MR-TOF systems without affecting MR-TOF resolution and mass accuracy.

Novel Aspect

High duty cycle multi-reflecting OA-TOF-MS using artefact-free multiplexed data acquisition method without compromising mass accuracy and resolving power.

References

1. A. Brock, N. Rodriguez, R. N. Zare, *Anal. Chem.*, 70, pp. 3735-3741, 1998.
2. J. Franzen, US Patent 2004164239, 26 Aug 2004.
3. M. Belov, C. A. Fancher, P. Foley, US Patent 6900431, 31 May 2005.
4. A. Verenchikov, WO Patent 135477, 3 Nov 2011.
5. P. M. Willis, J. Jaloszynski, US Patent 9786484, 10 Oct 2017.

TP-84 / SCREENING FOR ANTHRAQUINONES IN CHILEAN MUSHROOMS BY HPTLC - NEGATIVE ION DESI HIGH-RESOLUTION MASS SPECTROMETRY

Annegret Laub⁽¹⁾ - ***Ann-Katrin Sendatzki***⁽¹⁾ - ***Götz Palfner***⁽²⁾ - ***Norbert Arnold***⁽¹⁾ - ***Ludger Wessjohann***⁽¹⁾ - ***Jürgen Schmidt***⁽¹⁾

Leibniz Institute of Plant Biochemistry, Bioorganic Chemistry, Halle (1) - Universidad de Concepcion, Departamento de Botanica, Concepcion (2)

Keywords: HPTLC-negative ion DESI-HRMS/MS, mushrooms, anthraquinones

Introduction

Desorption electrospray ionization mass spectrometry (DESI-MS) is a powerful ambient ionization mass spectrometric technique [1]. Its coupling with high-performance thin-layer chromatography (HPTLC) provides a simple and robust methodological approach for separation and highly-sensitive detection of natural compounds in plants and fungi [2]. Moreover, high-resolution mass spectrometry (HRMS) ensures an outstanding selectivity and specificity of this technique. This methodology can be also used in natural product research for a fingerprint analysis in crude extracts [3]. We have applied a combination of a normal phase HPTLC with DESI-HRMS for profiling of anthraquinones in Chilean *Cortinarius* species.

Methods

Crude methanolic extracts of seven *Cortinarius* species were investigated by thin layer chromatography using 10x10 cm Silica 60 F254 HPTLC plates. The plates were developed in a horizontal CAMAG chamber using toluene/ethyl formate/ formic acid (10:5:3 v/v/v) as solvent system. The DESI-MS analysis was performed by an Orbitrap Elite mass spectrometer (Thermo Scientific, Bremen, Germany) combined with a 2D-DESI-source (Omnispray System OS-3201, Prosolia Inc., Indianapolis, USA) and operated in the negative ion mode. The MS and MS₂ spectra were acquired with the resolution of 30000 by scanning of individual TLC bands in the x-direction at a surface velocity of 200 µm/s. The data were evaluated using Xcalibur 2.2 SP1 software.

Results

We present HPTLC-negative ion DESI-HRMS results of the analysis of complex crude extracts from seven Chilean species of the genus *Cortinarius* concerning the occurrence of anthraquinones. These compounds are known as potential chemotaxonomic markers in fungi [4]. Using this method physcion, emodin, dermolutein, hypericin, skyrin, 7,7'-emodinphyscion, and flavomannin C were identified by their elemental composition and partly by their characteristic fragmentation patterns in comparison with reference compounds. It should be pointed out, that the HRMS approach also allows a mass spectral distinction of isobaric ions. This was demonstrated for the detection of some investigated anthraquinones whose nominal mass is accompanied by other compounds in the crude extracts.

Conclusions

HPTLC-negative ion DESI-HRMS was shown to be a powerful method for the analysis and profiling of anthraquinones in complex crude extracts of different Chilean fungi of the genus *Cortinarius*. Using this method six anthraquinones were identified by their elemental composition and partly by their characteristic fragmentation patterns in comparison with reference compounds. High-resolution mass spectrometry (HRMS) also allows a mass spectral distinction of isobaric ions which could be demonstrated for the detection of some investigated anthraquinones whose nominal mass is accompanied by other compounds in the crude extracts.

Novel Aspect

HPTLC-DESI-HRMS/MS as a powerful tool for chemotaxonomic analysis of Chilean mushrooms

References

1. Takats Z, Wiseman J M, Gologan B, Cooks R G, Science,306, 471-473 (2004).
2. Van Berkel G J, Tomkins B A, Kertesz V, Anal. Chem.,79, 2778-2789 (2007)
3. Bagatela B S, Lopes A P, Cabral C C, Perazzo F F, Ifa D R, Rapid Commun. Mass Spectrom.,29, 1530-1534 (2015)
4. Gill, M, Nat. Prod. Rep.,20, 615-639 (2003)

TP-85 / THE ADVANTAGES OF USING AN MRMHR ACQUISITION FOR INCREASED SENSITIVITY FOR TARGETED ANALYSIS USING A QTOF INSTRUMENT

Daniel McMillan⁽¹⁾ - ***Phil Taylor***⁽²⁾ – ***Katherine C. Hyland***⁽³⁾ - ***Chris Borton***⁽⁴⁾

SCIEX, Business Development, Warrington (1) - SCIEX, Marketing, Warrington (2) - SCIEX, Application Development, Redwood Shores (3) - SCIEX, Application Development, Framingham (4)

Keywords: QTOF, Targeted, MRM, Sensitivity

Introduction

When attempting to quantitate with a QTOF instrument using a fragment ion instead of the signal from the TOF acquisition is that the collision energy for a given fragment is not formed at an optimum collision energy. A novel feature of the SCIEX X500R QTOF System is the ability to acquire specific MRM pairs just as are acquired in a triple quadrupole instrument thus improving the sensitivity of the method. This presentation discusses the concept of an MRMHR experiment and demonstrates the advantage of this technique.

Methods

A set of standards were analyzed using a standard data dependent scan where the parent ion of each compound is selected for fragmentation in Q1 and the ion is fragmented using a collision energy of 35 volts. A quantitation method was created using the two most intense fragment ions and an instrumental

detection limit was determined. The same set of standards was re-analyzed using an MRMHR scan. The instrumental detection limit from the MRMHR data was determined and compared to the previous fragmentation data. The SCIEX X500R QTOF system has the unique ability to acquire MRM transitions with optimized collision energies resulting in triple quadrupole sensitivity but allows for the acquisition of full scan MS/MS for compound identification without sacrificing cycle time, and we will discuss the features of MRMHR

Conclusions

A limitation of traditional triple quadrupole mass spectrometers for quantitative applications is that compound identification relies on a ratio of two ions from the compound fragmentation. Interference from the matrix frequently causes this ratio to be unreliable for compound identification resulting in both false positive and false negative results. The slow speed of acquiring full scan MS/MS spectra from a triple quadrupole instrument prohibits the acquisition of useful spectra for compound identification due to degraded cycle times. One of the advantages of QTOF instruments is the speed of acquiring full scan MS/MS spectra during sample analysis without the slow cycle times observed for triple quadrupole instruments. Unfortunately, quantitation is accomplished using fragments from a full scan MS/MS where the collision energy is not optimized for a given transition.

Novel Aspect

This is the first demonstration of using a QTOF to acquire MRM transitions with optimized collision energies

TP-86 / CONTEXTUALIZATION OF DARK MATTER IN MS1 FT-ICR MS¹ DATA.

Franco Moritz^(1,2) - ***Moritz Kaling***⁽²⁾ - ***Jörg-Peter Schnitzler***⁽²⁾ - ***Philippe Schmitt-Kopplin***^(1,3)

Research Unit Analytical BioGeoChemistry, Department of Environmental Sciences, Helmholtz Zentrum München, IngolstaedterLandstrasse 1, 85764 Neuherberg, Germany (1) - Research Unit Environmental Simulation, Department of Environmental Sciences, Helmholtz Zentrum München, Ingolstaedter Landstrasse 1, 85764 Neuherberg, Germany (2) - Chair of Analytical Food Chemistry, Technische Universität München, Alte Akademie 10, 85354 Freising-Weihenstephan, Germany (3)

Keywords: FT-ICR-MS, Dark Matter, Annotation, Metabolomics

Introduction

Major proportions (70-95%) of ultra-high resolution MS-signals with known molecular formula cannot be found in metabolomic databases. This mass spectrometric dark matter [1, 2, 3] exhibits statistical significance and is usually not identifiable via MS/MS. This lack of biological MS1-signal contextualization deems major parts of data useless. Interpretation methods that go beyond mere m/z -to-database matching are required to tackle this problem.

Methods

Ultra-high resolution MS data on *Populus X. canescens* were annotated by virtue of mass difference network propagation. Mass difference enrichment analyses (MDEAs), mapping to databases via MDs and compositional similarity and PLS regression were used to extract molecular patterns in association to Isoprene synthase knockdown and experimental perturbations. SIM-STITCH MS2 scans were employed to investigate the plausibility of predictions [4].

Results

4,335 m/z peaks that obtained molecular formula annotation occurred non-randomly across 211 mass spectra. Merely 129 molecular formulas (3%) matched to KEGG pathways, while several hundreds of peaks discriminated wild type plants from ISPS knock-downs. MDEAs highlighted the molecular building blocks

that were used to build up either class's biomass. The usage patterns of the building blocks implied by MS1 were consistent with SIM-STITCH MS2 results. MDEA confirmed prior knowledge on the effects of ISPS knock-downs. The highlighted mass difference building blocks (MDBs) were used to connect the mass spectrometric dark matter to KEGG database entries, which resulted in consistent transformation patterns, explaining how the different genotypes transform KEGG compound classes and metabolites from other databases into biomass. Compound class conserving MDs were elaborated. Compound class distributions were in-line with results from compositional similarity approaches with the power to contextualize almost 100% of the obtained data.

Conclusions

The mass-difference space in MS1 spectra contains valuable information on usage patterns of molecular building blocks. MDEA can highlight the processes that drive the synthesis of biomass after fundamental perturbations to systemic carbon flux. Dark matter contextualization will help narrowing down fruitful fields of research by highlighting compound classes that change systematically.

Novel Aspect

UHR-MS peaks are considered as being altered versions of well-known metabolites. We introduce parametric- and non-parametric ways to search for such systemic alterations.

References

1. da Silva RR, Dorrestein PC, Quinn RA. Illuminating the dark matter in metabolomics. *Proc Natl Acad Sci USA* (2015). 112(41), 12540-12550.
2. Joseph B, Corwin JA, Baohua L, Atwell S, Kliebenstein DJ. Cytoplasmic genetic variation and extensive cytonuclear interactions influence natural variation in the metabolome. (2013) *eLife* 2e00776.
3. Dührkop K, Shen H, Meusel M, Rousu J, Böcker S. Searching molecular structuredatabases with tandem mass spectra using CSI:Finger ID. (2015) *Proc Natl Acad Sci USA* 112(41) 12580–12585.
4. Moritz F, Kaling M, Schnitzler J-P, Schmitt-Kopplin Ph. Characterization of poplar metabolotypes via mass difference enrichment analysis. (2017) *Plant, Cell & Environment* (2017) 40, 1057-1073.

TP-87 / SIMULTANEOUS DETERMINATION OF NINE WATER-SOLUBLE SYNTHETIC COLORANTS IN COSMETICS BY LC-Q-TOF/MS AND LC-MS/MS

Seongsoo Park - Jun Hyung Kim - Gihaeng Kang - Ayoung Min - Ji Hyun Lee - Hyoung-Joon Park - Taeik Hwang - Seok Heo - Kwon Yong Eom - Dong Woo Shin - Sung-Kwan Park - Hoil Kang

National Institute of Food and Drug Safety Evaluation, Division of Advanced Analysis, Chungcheongbuk-do

Keywords: synthetic colorants, cosmetic, LC-Q-TOF/MS, LC-MS/MS

Introduction

In this study, LC-Q-TOF/MS and LC-MS/MS method were developed to identify and measure the illegally adulterated synthetic colorants (Acid yellow 9, Brilliant black BN, Orange GGN, Yellow 2G, Red 2G, Green S, Carmoisine, Erythrosin B, Patent blue V) that may be found in cosmetics circulated in Korea. .

Methods

We have established and optimized the simultaneous analysis of nine water-soluble synthetic colorants in facial toner using LC-Q-TOF/MS and LC-MS/MS. Colorants-spiked facial toner samples (0.2 g) were prepared with 50 mL of 20% methanol:chloroform (9:1) into a Teflon centrifuge tube. Sample extraction was performed by vortexing the sample for 2 min and ultrasonication for 30 min. The extract was centrifuged for 5 min at 3000 rpm and upper layer was collected in volumetric flask of 50 mL. All solutions were filtered through a 0.22 µm syringe filter (PVDF) prior to LC-Q-TOF/MS and LC-MS/MS analysis.

Results

The method was validated for specificity, linearity, limit of detection (LOD), limit of quantification (LOQ) and recovery. The linear correlation coefficients were good (>0.99), the limits of detection and quantitation of the method ranged from 0.0039 to 0.0625 $\mu\text{g/mL}$ and from 0.0117 to 0.1875 $\mu\text{g/mL}$, respectively. The recovery for this LC-MS/MS procedure was 82.5–104.7%, with an RSD less than 15% ($n=3$).

Conclusions

In this study, a fast, simple, and accurate method was developed for the identification and quantification of nine water-soluble synthetic colorants in facial toner by using LC-Q-TOF/MS and LC-MS/MS. The method was completely validated showing satisfactory data for specificity, linearity, limit of detection, limit of quantification and recover. Therefore, this study will be helpful for comprehensive analysis and monitoring synthetic colorants in cosmetics.

Novel Aspect

In this study, we have developed analytical methods for illegally adulterated synthetic colorant in cosmetics circulated in Korea. The high sensitivity, high resolution and accurate mass determination, LC-Q-TOF/MS method is a useful analytical tool for synthetic colorants. Also, the results indicated that the method can be used for rapid and accurate screening of synthetic colorants present in cosmetics

References

1. Eugenia G., Maria C., Pablo L., Maria L., Carmen G., Journal of Chromatography A, 1415, 27-37 (2015)
2. Yanping X., Yuluan W., Xindong G., Yujing K., Halying L., Donghui L., Yiguanq C., Anal. Methods, 5, 1965 (2013)

TP-88 / SFC-ESI-QTOF-MS: DEVELOPING AN INNOVATIVE ONE-RUN METHOD FOR CHOCOLATE QUALITY AND AUTHENTICITY CONTROL

Michaela Rektorisova - Vojtech Hrbek - Monika Tomaniova - Jana Hajslova

University of Chemistry and Technology, Prague, Department of Food Analysis and Nutrition, Prague

Keywords: High-resolution mass spectrometry, Parallel reaction monitoring, Triacylglycerols, Cocoa butter equivalents, Cocoa solids

Introduction

To assess chocolate quality and authenticity, a combination of various analytical procedures is involved, thereby making the process time-consuming and costly. [1] Thus, we investigated the potential of supercritical fluid chromatography coupled to quadrupole-time of flight mass spectrometry (SFC-QTOF-MS) to replace 'classic' methods. Its specific separation character and selective detection make this technique suitable for authentication of fatty foods.

Methods

Samples were prepared by sequential extraction using hexane in the first step, followed by water extraction of the residue. Analysis of combined extracts was performed by Waters SFC-QTOF-MS with electrospray ionization (ESI). To obtain most comprehensive information about sample composition, the MS method comprised a full MS1 scan to detect potential unpredicted compounds, and several time-scheduled targeted MS/MS functions (parallel reaction monitoring).

Results

A single 8 minute analytical run enabled us (i) to determine the presence of cocoa butter equivalents and milk fat according to triacylglycerols, (ii) to calculate dry non-fat cocoa solids based on theobromine and caffeine content, and (iii) to screen sugars, *i.e.* the key components involved in legislation. Optimizing

parameters in the individual MS/MS functions (scan time, target enhancement, collision energy) led to a significant adjustment in analyte responses according to their concentration range in samples, *e.g.* decreased response of 3 major triacylglycerols of cocoa butter in contrast to increased response of minor triacylglycerol specific for milk fat. Attention had to be paid to background contamination (*e.g.* phthalates) to avoid the suppression of target ion intensity. Thus, in addition to an effort to minimize such contamination, it was useful to acquire a full MS scan not only to obtain comprehensive information about samples, but also to prevent quantification errors caused by abrupt changes in background intensity.

Conclusions

SFC-QTOF-MS has proved to be a powerful tool to assess chocolate quality and authenticity, enabling fast, simultaneous determination of cocoa butter equivalents, non-fat cocoa solids and sugars. We have demonstrated a unique separation (compared to reversed-phase liquid chromatography) and advantages of using parallel reaction monitoring for targeted analysis. In addition, we have shown the importance of acquiring a full MS scan to aid interpreting the results.

Novel Aspect

Fast and cost-effective method for the assessment of several chocolate quality and authenticity parameters that are usually analyzed by a combination of standard methods.

References

1. Joint FAO/WHO Codex Alimentarius Commission. Codex Alimentarius, Codex Stan 87-1981(2016).

TP-89 / USING ATOMIC LAYER DEPOSITION (ALD) TO ENHANCE MCPS FOR MASS SPECTROMETRY

Matthew Breuer⁽¹⁾ - **Paula Holmes**⁽¹⁾ - **John Harper**⁽²⁾

Photonis, Research and Development, Sturbridge, MA (1) - Photonis, Research and Development, Brive (2)

Keywords: Microchannel plate (MCP), Atomic Layer Deposition (ALD), detector

Introduction

A microchannel plate (MCP) is an array of identical parallel open channels with diameters typically between 5 and 25 microns in a glass plate with a thickness that is typically 40-60 times the channel diameter. MCPs are made of a special lead-silicate glass that, when properly processed, makes the MCP an excellent ion detector and electron multiplier, capable of high gain and low noise. While it has long been possible to coat the input and output surfaces of MCPs to enhance their performance, historically it has been prohibitively difficult to use a coating to affect the inner channel surface that is responsible for electron multiplication. Atomic layer deposition (ALD), which is a process that builds conformal coatings one atomic layer at a time using a sequence of self-terminating gas-solid reactions, offers an opportunity to coat these high aspect ratio surfaces with high uniformity and enhance some aspects of the MCP performance as a result.

Methods

Standard and coated MCPs were fabricated and their relevant performance characteristics were characterized using industry standard test methods for MCPs. These characteristics include overall gain, gain stability (lifetime), and linear current output. In addition, coated and uncoated MCPs were assembled in detector assemblies intended for mass spectrometry and their instrument-specific characteristics were evaluated using the instrument tune sequences.

Results

Coatings of high electron yield materials such as aluminum oxide can result in higher overall gains when the MCPs are operated at high voltage and lower gains when the MCPs are operated at lower voltages. This is

in line with expectations based on the dependence of the kinetic secondary electron emission yield curves as a function of primary electron energy. Coated MCPs have good gain stability during operation at typical vacuum levels (<1E-5 Torr) similar to uncoated MCPs. When used in time of flight detector assemblies, aluminum oxide coated MCPs generally tune at a lower voltage. Aluminum oxide coatings can also be used to form an inert barrier layer on standard MCPs. This can affect long-term resistance stability when the MCPs are exposed to air or dry nitrogen. It is typical for the overall resistance of uncoated MCPs to increase gradually over time. Coating of the MCPs can reduce this effect.

Conclusions

ALD can be used to uniformly coat the high-aspect ratio inner channel surfaces of MCPs and affect their electron multiplication behavior. This provides an additional method for tailoring MCPs to meet the required performance characteristics of particular applications.

Novel Aspect

An alternate electron multiplication coating can be applied to MCPs to enhance performance.

TP-90 / THE COMBINATION OF HYDROPHILIC INTERACTION LIQUID CHROMATOGRAPHY AND ION MOBILITY MASS SPECTROMETRY FOR DIESEL FUEL ANALYSIS

Elize Smit⁽¹⁾ - *Stefan De Goede*⁽²⁾ - *Egmont Rohwer*⁽³⁾

University of Johannesburg, Department of Chemistry, Johannesburg (1) - SASOL Energy, Energy Technology, Sasolburg (2) - University of Pretoria, Department of Chemistry, Pretoria (3)

Keywords: High resolution mass spectrometry, Hyphenated techniques, Hydrophilic Interaction Liquid Chromatography, Ion Mobility, Petroleomics

Introduction

The characterization of trace polar species in fuels is important due to their influence on the fuel's physical properties, which has given rise to the field of petroleomics [1]. Due to the complex nature of diesel and other crude-derived fuels its chemical composition is often investigated by considering classes of compounds rather than individual chemical species [2,3,4].

Methods

The chromatographic separation of diesel fuel extracts using hydrophilic interaction liquid chromatography (HILIC) was combined with high resolution mass spectrometry (TOFMS). Selective ionization of trace polar species was achieved with electrospray ionization and atmospheric pressure chemical ionization. Furthermore, ion mobility mass spectrometry was used to obtain additional chemical information.

Results

The separation of the polar fraction of diesel into classes of compounds, *i.e.* based on functionality, was achieved. Complementary information was obtained by considering both positive and negative ion mass spectra. Separation of species that readily formed positive ions was influenced by heteroatom content as well as degree of unsaturation (or DBE). Visualisation of the multidimensional data was explored in order to facilitate the interpretation of chemical information obtained. More than half of the detected m/z values were observed to show isomers separated by drift and/or retention time. A large amount of structural information is therefore gained by pre-separation despite the loss in sensitivity. A novel graphical representation has been introduced where Kendrick mass defect (instead of m/z values) is plotted against drift time to visualise information regarding heteroatom content and double bond equivalents (DBE) against structural information, resulting in an orthogonal system.

Conclusions

Hydrophilic interaction liquid chromatography was used to separate methanol extracts of diesel fuel in terms of heteroatom content (*i.e.* group-type separation). Ion mobility mass spectrometry provided additional chemical information of species in diesel fuel extracts. The multidimensional data obtained from hyphenated techniques can be visualized in unique ways to facilitate interpretations.

Novel Aspect

The group-type separation of diesel fuel extracts was achieved with hydrophilic interaction chromatography. Mobilograms were modified to incorporate mass defects.

References

1. Marshall A. G., Rodgers R. P., *Acc. Chem. Res.* 37, 53–59 (2004).
2. Blomberg J., Schoenmakers P. J., Beens J., Tijssen R., *J. High Resolut. Chromatogr.* 20, 539–544 (1997).
3. Dallüge J., Beens J., Brinkman U. A. T., *J. Chromatogr. A* 1000, 69–108 (2003).
4. Hughey C. A., Hendrickson C. L., Rodgers R. P., Marshall A. G., Qian K., *Anal. Chem.* 73, 4676–4681 (2001).

TP-91 / APPLYING HIGH-RESOLUTION MASS SPECTROMETRY TO ELUCIDATE THE FORMATION AND FATE OF OZONATION TPS IN WASTEWATER TREATMENT

Jennifer Schollee⁽¹⁾ - *Rebekka Gulde*⁽¹⁾ - *Marc Bourgin*⁽¹⁾ - *Moreno Rutsch*⁽¹⁾ - *Julian Fleiner*⁽²⁾ - *Marc Boehler*⁽²⁾ - *Juliane Hollender*⁽¹⁾ - *Urs von Gunten*⁽³⁾ - *Christa McArdell*⁽¹⁾

Eawag, Environmental Chemistry, Duebendorf (1) - Eawag, Process Engineering, Duebendorf (2) - Eawag, Water Resources & Drinking Water, Duebendorf (3)

Keywords: Non-target screening, Ozonation, Transformation products, Advanced wastewater treatment

Introduction

Ozonation has increasingly been applied for micropollutant abatement in wastewater and high removal rates have been observed for most known micropollutants. However, the formation and fate of ozonation transformation products (OTPs) is not well understood. High resolution tandem mass spectrometry (HRMS/MS) with non-target analysis was applied to study unknown OTPs in real wastewaters in Switzerland.

Methods

Samples were collected from two wastewater treatment plants (WWTPs) before and after ozonation and measured with Orbitrap HRMS/MS coupled to reverse-phase liquid chromatography (RPLC). Hierarchical cluster analysis grouped features based on intensity patterns and the fate of likely OTPs was investigated. Additionally, a method was developed for HRMS coupled to hydrophilic interaction liquid chromatography (HILIC) to measure highly polar OTPs.

Results

Out of approximately 8000 HRMS features detected at one WWTP, 600-700 features (9%) were determined to be likely OTPs based on their intensity patterns. In subsequent post-treatment, 54-83% of likely OTP features were well-removed, with granular activated carbon (GAC) filters (run up to 40,000 bed volumes) performing slightly better than the biofilters investigated. It is speculated that sorption in the GAC filters accounted for the higher OTP abatement compared to biofilters, where biodegradation was the primary mechanism for OTP abatement. It was also found that OTPs that were well-removed in post-treatment had significantly higher retention times in RPLC ($p < 0.05$, Mann-Whitney U-test), suggesting recalcitrance of polar OTPs. To get information on highly polar OTPs, a new method combining HILIC with HRMS/MS was developed and preliminary results indicate that previously undetected OTPs can now be investigated with this method. Formation of OTPs from 50 micropollutants is being further studied with HILIC- and RPLC-HRMS/MS.

Conclusions

The combination of HRMS with non-target analysis provides a more comprehensive picture of the formation and abatement of OTPs. By developing a HILIC-HRMS/MS method, unknown highly polar OTPs can also be evaluated. The ability to monitor OTPs in WWTPs is critical to elucidate persistent OTPs that might make their way to the aquatic environment.

Novel Aspect

Novel analytical and data analysis methods applied to HRMS/MS data were used to evaluate the formation and abatement of OTPs during advanced wastewater treatment.

TP-92 / ANALYSIS OF CATIONIC POLAR PESTICIDES IN BABY FOOD USING ION CHROMATOGRAPHY AND MS/MS OR HIGH RESOLUTION ACCURATE MASS SPECTROMETRY

Frans Schoutsen⁽¹⁾ - *John E Madden*⁽²⁾ - *Jeffrey S Rohrer*⁽²⁾ - *Jonathan Beck*⁽³⁾

Thermo, Takkebijsters 1, Breda (1) - TFS, CA, Sunnyvale (2) - TFS, CA, San Jose (3)

Keywords: LCMS, HRAM, Ion Chromatography, Chlormequat, QuPPE

Introduction

In this publication we demonstrate analysis of quaternary amine pesticides in baby food samples using cation-exchange chromatography. Polar cationic pesticides such as chlormequat, diquat, mepiquat, paraquat, morpholine, trimethylsulfonium, and propamocarb may occur as residues in food, but may be excluded in pesticide monitoring programs due to challenges using generic multi-residue methods.

Methods

Cation-exchange chromatography with serial detection by suppressed conductivity and mass spectrometry. These techniques are demonstrated by targeted MS/MS on a triple quadrupole mass spectrometer and by high resolution accurate mass spectrometry (HRAM MS)

Results

A chromatographic separation was demonstrated on a prototype column. The six pesticides had good accurate mass, meeting the SANTE requirements of <5 ppm. Native cationic pesticides of interest were not found in the samples. Sensitivities in the low µg/L were determined by spiking pesticides in the samples. Good accuracy was found, with recoveries between 80-120%.

Conclusions

Ion chromatography coupled with HRAM mass spectrometry demonstrated clear benefits for challenging cationic pesticide separation and identification.

Novel Aspect

Ion chromatography coupled with HRAM mass spectrometry demonstrated clear benefits for challenging cationic pesticide separation and identification.

TP-93 / AN AUTOMATED WORKFLOW FOR HIGH-THROUGHPUT NICOTINE QUANTITATION USING HIGH PERFORMANCE LIQUID CHROMATOGRAPHY COUPLED TO HIGH RESOLUTION ACCURATE MASS SPECTROMETRY FOR PROTOTYPE TESTING

Sandra Sendyk - *Quentin Dutertre* - *Arno Knorr* - *Julia Hoeng* - *Mark Bentley*

Philip Morris, Science & Innovation, Neuchatel

Keywords: Mass spectrometry, high-throughput, automation, method validation

Introduction

Human organotypical 3D tissue cultures are routinely used for in vitro toxicological assessment of the biological impact upon exposure to smoke from a reference cigarette 3R4F and aerosol from heat-not-burn products at the air-liquid interface using the Vitrocell® system. Nicotine, one of the major constituent of tobacco can be used as biomarker of exposure.

Methods

Cigarette smoke or aerosol of heat-not-burn products was deposited at Phosphate Buffer Saline solution (PBS), serving as a surrogate matrix in the Vitrocell® system. A reverse phase high-performance liquid chromatography coupled to high resolution accurate mass spectrometry was used.

The procedure was automated using a protocol developed in Pipeline Pilot (Biovia®) from acquisition to reporting.

Results

The method was designed enabling analysis of a large number of samples within a short time.

Selectivity of chromatographic separation was proven for nicotine and its internal standard (nicotine-methyl-d₃), separated from other interferences.

A high-throughput method was designed and enables quantitation of nicotine within a dynamic range of three order of magnitude. The method was validated within a working range from 50 ng/mL to 2000 ng/mL in PBS. Linearity of the method was achieved by a 8-point calibration curve. Furthermore accuracy, precision and instrumental repeatability were assessed.

Quantitation was performed by isotopic dilution technique using stable isotope labelled nicotine as internal standard. The method is reliable for the analysis of subsequent sample batches and demonstrated high robustness in practical applications.

Conclusions

An automated workflow for high-throughput nicotine determination was established delivering robust and accurate data. The method was validated testing accuracy, linearity, instrument repeatability and stability matching acceptance criterion according to standard regulation [1, 2]. It shows suitable performances for routine based application to accurately determine nicotine concentration in PBS.

Novel Aspect

An automated workflow for high-throughput PBS exposed samples analysis by liquid chromatography coupled to high resolution accurate mass spectrometry.

References

1. Committee for Medicinal Products for Human Use, European Medicines Agency Science Medicines Health, 2011
2. Kalakuntla Rao R., Santosh Kumar K., Journal of Pharmaceutical Sciences and Research, 1, 1-10, 2009

TP-94 / CHARACTERIZATION OF STRESS DEGRADATION PRODUCTS OF FOSAMPRENAVIR CALCIUM USING Q EXACTIVE™ PLUS ORBITRAP AND LTQ XL™ MASS SPECTROMETERS, AND PREDICTION OF THEIR TOXICITY USING ADMET PREDICTOR™

Dilip Kumar Singh⁽¹⁾ - ***Mayurbhai Kathadbhai Ladumor***⁽¹⁾ - ***Chinmaya Narayana Kotimoole***⁽²⁾ - ***Kedar Balaji Batkulwar***⁽²⁾ - ***Abhijeet Yashwantrao Deshpande***⁽²⁾ - ***Sanjeev Giri***⁽²⁾ - ***Saranjit Singh***⁽¹⁾

Department of Pharmaceutical Analysis, National Institute of Pharmaceutical Education and Research (NIPER), S. A. S. NAGAR (1) - DMPK, Aurigene Discovery Technologies Limited, Hyderabad (2)

Keywords: Fosamprenavir calcium, Stress studies, Degradation products, Orbitrap, ADMET Predictor™

Introduction

Stress testing is an essential part of a regulatory-compliant stability program for all types of drugs [1]. It involves i) subjecting the drug to a variety of conditions, ii) separation of degradation products (DPs) by a stability-indicating method (SIM), and iii) characterization of DPs with the help of advanced spectroscopy tools [2]. These studies were systematically accomplished on fosamprenavir calcium (FPV), a protease inhibitor used to cure HIV [3].

Methods

The drug was subjected to stress under hydrolytic, oxidative, photolytic and thermal conditions. It proved stable in the latter three conditions. The hydrolytic DPs were separated using a SIM, and were initially characterized from mass results obtained by employing UHPLC-Thermo Q-Exactive™ Plus Orbitrap and LC-LTQ XL™ spectrometers. Toxicity of each DP was predicted using ADMET Predictor™ software (Simulation Plus).

Results

Six DPs (DP 1-DP 6) were formed on acid and base stress. DP 1, DP 2, and DP 4-DP 6 were generated in acid (2N HCl, 75°C, 6 d), while DP 3 was the lone DP in basic (2N NaOH, 75°C, 1 d) condition. DP 2 was a common product. Mass studies revealed that the drug was protonated at two different sites, each following a separate fragmentation route. HRMS and MS_n data allowed suggesting comprehensive fragmentation pathways of the drug as well as its DPs. The information allowed successful delineation of structures of all DPs. Major hydrolytic routes were cleavage of carbamate ester of the drug to form DP 2, and dephosphorylation of both the drug and DP 2 to form DP 5 and DP 1, respectively. Methylation of the drug in presence of methanol resulted in DP 6 (at aniline nitrogen) and DP 3 (upon replacement of THF moiety). DP 6 underwent carbamate ester cleavage to result in DP 4. Toxicity predictions highlighted that DP 1 was associated with propensity for phospholipidosis, DP 2-DP 4 had susceptibility towards skin sensitization, while DP 5 was associated with reproductive toxicity.

Conclusions

Stress degradation study on FPV provided significant information on its degradation behavior. The drug was found to be susceptible only to acid/base hydrolysis. A total of six DPs were formed, which were characterized using HRMS and MS_n data. The DPs were also evaluated for their in silico toxicity using ADMET Predictor™ software.

Novel Aspect

Degradation profile of FPV was established and justified through relevant mechanistic explanations. In silico prediction studies highlighted added toxicities of DPs.

References

1. Singh S., Junwal M., Modhe G., Tiwari H., Kurmi M., Parashar N., Sidduri P., TrAC Trends in Analytical Chemistry, 49, 71-88 (2013).
2. Bakshi M., Singh S., Journal of Pharmaceutical and Biomedical Analysis, 28, 1011-1040 (2002).
3. Pokorná J., Machala L., Řezáčová P., Konvalinka J., Viruses, 1, 1209-1239 (2009).

TP-95 / QUADRUPOLE-RESOLVED ALL IONS MS/MS FOR RELIABLE QUANTITATION OF PESTICIDES IN COMPLEX MATRICES USING UNTARGETED ACQUISITION

László Tölgyesi⁽¹⁾ - Andreas Mielcarek⁽¹⁾ - Christian Klein⁽²⁾ - Tom Knotts⁽²⁾ - Dorothy Yang⁽²⁾ - Bill Frazer⁽²⁾ - Tarun Anumol⁽³⁾ - John Lee⁽⁴⁾ - William E. Barry⁽²⁾

Agilent Technologies Sales & Services GmbH & Co. KG, Global Laboratory Solution Sales (LSS) Marketing, Waldbronn (1) - Agilent Technologies, Inc., R&D, Santa Clara, CA (2) - Agilent Technologies, Inc., Segment Marketing, Little Falls, DE (3) - Agilent Technologies, Inc., Segment Marketing, Bracknell, Berkshire (4)

Keywords: QTOF, untargeted screening, pesticides, HRMS

Introduction

For pesticides screening, untargeted data collection utilizing LC/MS High Resolution Mass Spectrometry allows labs to continually evolve the scope of targets without the need to change a base method. Measurement of fragment ions can sometimes suffer from interferences, effecting the reliability of ion ratio and ultimately the ability to detect much reduced fragments. To overcome this, precursors are resolved into smaller bands of ions on a prototype QTOF.

Methods

The quadrupole isolation windows are adjustable by changing DAC values ranging from 0 to 255, effecting the applied RF values. By using calibrant ions, the isolation window width and rectangularity could be determined. Multiple acquisition methods were generated, dividing the mass range into sections between 10 amu and up to several hundred amu. For proof of concept, avocado and black tea QuEChERS extracts, spiked to 2 and 10 ppb with pesticides, were evaluated.

Results

The rectangularity of the quadrupole window width needed to be determined, because the shape of the isolated window determines the necessary overlap when moving from one window to the next. Generating empirically the needed DAC values for a 30 amu and a 50 amu window, and scanning the quad in 0.1Da steps, less than 1 amu was needed to reach FWHM on both sides of the isolation window. Calibrant mass 322 was used to establish beginning, end and half maximum of the window.

In the case of Imazalil, however, the All Ion(AI) spectra showed a matrix interference of one of its fragment ions, leading to a mass shift as well as an EIC height which was overlaid with the matrix contribution. The effect of matrix contribution went up with lower concentration of the pesticide. By using 30 amu isolation window width, the matrix interference was completely removed, allowing the reliable quantitation of Imazalil in both avocado and black tea matrices with identical ion ratio compared to the standard in neat solvent matrix.

Conclusions

Comparing the analyses of AI MS/MS alone with quadrupole resolved All Ions MS/MS on the spiked pesticide samples, we have noticed that the fragmentation spectra looked substantially cleaner compared to All Ions MS/MS spectra. Nevertheless, in multiple cases ion-ratio of fragment ions remained the same. The dependence of pesticide concentration, quadrupole window width (fixed and variable) as well as matrix complicity will be further evaluated and presented.

Novel Aspect

Increased confidence in compound qualification and quantitation on a Q-TOF mass spectrometer by utilizing TQ-like post-processing criteria using untargeted acquisition.

TP-96 / IDENTIFICATION OF ELECTROCHEMICAL DEGRADATION PRODUCTS OF SELECTED PHARMACEUTICALS BY HIGH RESOLUTION MASS SPECTROMETRY

Andrea Vojs Staňová^(1,2) - ***Erika Medvecka***⁽¹⁾ - ***Adam Borik***⁽²⁾ - ***Roman Grabic***⁽²⁾ - ***Pavel Škvára***⁽¹⁾ - ***Tomas Mackulak***⁽³⁾ - ***Marian Marton***⁽⁴⁾ - ***Marian Vojs***⁽⁴⁾

Comenius University in Bratislava, Faculty of Natural Sciences, Department of Analytical Chemistry,

Bratislava, Slovak Republic (1) - University of South Bohemia in České Budějovice, Faculty of Fisheries and Protection of Waters, South Bohemian Research Center of Aquaculture and Biodiversity of Hydrocenoses, Vodňany, Czech Republic (2) - Slovak University of Technology in Bratislava, Faculty of Chemical and Food Technology, Department of Environmental Engineering, Bratislava, Slovakia (3) - Slovak University of Technology in Bratislava, Faculty of Electrical Engineering and Information Technology, Institute of Electronics and Photonics, Bratislava, Slovak Republic (4)

Keywords: mass spectrometry, boron doped diamond electrodes, pharmaceuticals, wastewater treatment

Introduction

In recent years, pharmaceuticals are becoming more frequently found in the aquatic environment [1] which indicates limitations of their removal in conventional wastewater treatment plants [2]. Therefore, new treatment processes that are more efficient and not harmful for the environment are developed. Among different technologies for advanced water treatment belongs boron doped diamond (BDD) [3,4].

Methods

Residual pharmaceuticals and their metabolites produced during the electrochemical degradation process were determined using TSQ Quantiva LC/MS/MS (Thermo Fisher Scientific, San Jose, CA, USA). For the identification of electrochemical degradation products high-resolution mass spectrometers Q-Exactive MS (Thermo Fisher Scientific) and IT-TOF-MS (Shimadzu, Kyoto, Japan) were used.

Results

Electrochemical degradation of selected pharmaceuticals was performed in electrochemical cell with two parallel BDD electrodes. Based on our preliminary results, different types of BDD electrodes (grown with different deposition parameters, like CH_4/H_2 , B/C in gas phase) were used for water treatment in the first part of this study. Obtained results indicate that this method is able to effectively degrade a wide spectrum of pharmaceuticals that are insufficiently removed by traditional processes at WWTP. Micropollutants such as caffeine, sulfamethoxazole, diclofenac, carbamazepine, valsartan, trimethoprim and venlafaxine were degraded with more than 50% efficiency. In the second part of our study we focused on the identification of electrochemical degradation products. To identify potential electrochemical degradation products two different high resolution mass spectrometers were used.

Conclusions

This work deals with the possibilities of boron doped diamond electrodes for electrochemical degradation of different pharmaceuticals in wastewater. Using high resolution mass spectrometry, it was possible to identify the products of electrochemical degradation of selected pharmaceuticals.

Acknowledgements

This work was supported by the grant of the Slovak Research and Development Agency (APVV-16-0124).

Novel Aspect

Identification of electrochemical degradation products of selected pharmaceuticals by high resolution mass spectrometry.

References

1. Grabic R., Fick J., Lindberg R. H., Fedorova G., Tysklind M., Talanta, 100, 183-195 (2012).
2. Golovko O., Kumar V., Fedorova G., Randak T., Grabic R., Environ. Sci. Pollut R, 21, 7578-7585 (2014)
3. Mackulák T., Marton M., Radičová M., Vojs Staňová A., Grabic R., Bírošová L., Nagyová K., Vojs M., Bodík I., Brandeburová P., Gál M., Monats. Chem., 148, 539-548 (2017).
4. Rajab M., Heim C., Letzel T., Drewes J. E., Helmreich B., Int. J. Environ. Sci. Technol. 13, 2539-2548 (2016).

TP-97 / REMOVAL OF ISOBARIC INTERFERENCES BY LASER RESONANCE IONIZATION MASS SPECTROMETRY FOR DETERMINATION OF $^{126}\text{Sn}/^{121\text{m}}\text{Sn}$ RATIO

Wenliang Wang - Xiaopan Shen - Lihua Zhai - Hu Deng - Zhiming Li

Northwest Institute of Nuclear Technology, applied physics and chemistry department, Xi'an

Keywords: Laser resonance ionization mass spectrometry; Isobaric interference; Isotope ratio; Fission product; Spent fuel reprocessing solution

Introduction

As fission product, ^{126}Sn and $^{121\text{m}}\text{Sn}$ are considered as important nuclides to monitor the nuclear contamination [0-0]. The common methods for determination of $^{126}\text{Sn}/^{121\text{m}}\text{Sn}$ ratio are inductively coupled plasma mass spectrometry (ICP-MS) [0], thermal ionization mass spectrometry (TIMS)[0], and accelerator mass spectrometry (AMS)[0]. However, the isobaric interferences caused by ^{126}Te and ^{121}Sb are the main difficult problem in the those methods. A complicated sample preparation chemistry procedure is commonly adopted to separate the interfering species, which may cause waste of sample and introduce new interferences.

Laser resonance ionization mass spectrometry (LRIMS) has both advantages of laser resonance ionization spectroscopy and mass spectrometry. By selectively ionizing only one element at a time through a resonance excitation process, isobaric interferences are essentially eliminated [0]. In this paper, a LRIMS method of measurement of $^{126}\text{Sn}/^{121\text{m}}\text{Sn}$ ratio in the sample extracted from the spent fuel reprocessing solution was proposed. The removal of isobaric interferences in this method was confirmed by analysis of a tin-antimony-tellurium mixture sample.

Methods

Our LRIMS device developed indigenously mainly includes laser system, atomization source, electromagnetic analyzer, data acquisition and processing system. The basic process of the proposed LRIMS method is as follows: firstly, the sample is atomized by electric heating; secondly, the atoms are selectively ionized by the dye lasers with specific wavelengths; lastly, the ions are draw into the electromagnetic analyzer for isotopic ratio analysis.

The laser system consists of three tunable dye lasers pumped by 532 nm Nd:YAG solid-state lasers with pulse frequency of 10 kHz. The tunable wavelength range is 550 ~ 850 nm, which can be extended to shorter wavelength (275 to 425 nm) by frequency doubling method.

The atomization source is designed by electric heating tantalum filament with 18 mm long, 0.7 mm width, and 0.045 mm thickness . A grid electrode is mounted behind the atomization source to shield the ions produced by the thermal surface ionization. The detection efficiency of the sample with and without grapheme oxide solution are compared.

Three laser beams are coupled into the 3×1 fiber combiner, and focused by achromatic lens. Through carefully adjustment, the special overlap between the laser beam and atomic beam is maximized. The electromagnetic analyzer with the mass resolution of 500 and the abundance sensitivity of 2×10^{-6} is used to measure tin isotopic ratios.

In order to evaluate the isobaric interference in LRIMS method, a mixed sample of Sn, Sb and Te with a ratio of 1:1:1 is prepared to evaluate the isobaric interference. The $^{126}\text{Sn}/^{121\text{m}}\text{Sn}$ isotopic ratio in a tin sample extracted from the spent fuel reprocessing solution is determined.

Results

The auto-ionization energy spectrum of tin atom is obtained by laser resonance ionization spectrometry. According to the LRIMS signal of tin ions as a function of each laser wavelength, an efficient three-color-three-photon resonance ionization scheme of tin is confirmed. The laser wavelengths of each excitation/ionization step are optimized as $\lambda_1=286.42$ nm, $\lambda_2=811.63$ nm, $\lambda_3=823.72$ nm. In order to achieve saturated resonance excitation and ionization, the laser average power are set as 18 mW, 230 mW,

and 1700 mW respectively. The linewidth of each laser is tuned as about 10 pm to reduce the isotopic discrimination effect.

The temperature of the filament required for adequate neutral tin species production is optimized as 600-700°C. The total detection efficiency of 1 µg tin sample is measured as about 3×10^{-5} by mixing samples with grapheme oxide solution, which is as 4.5 times higher as the sampling without grapheme oxide solution. The results show that a layer of grapheme oxide gives better control of evaporation and reduces loss of sample as molecular oxide species.

The tin isotopic ratios in a tin-antimony-tellurium mixture with a ratio of 1:1:1 are measured. The major tin isotope ratios in the sample are determined and the relative standard deviations of $^{116}\text{Sn}/^{120}\text{Sn}$, $^{117}\text{Sn}/^{120}\text{Sn}$, $^{118}\text{Sn}/^{120}\text{Sn}$, $^{119}\text{Sn}/^{120}\text{Sn}$ are all less than 1%. There are no signal at $m/z=121$ (^{121}Sb) and 126 (^{126}Te) during the experiment indicating a selectivity factor for Sn exceeding 108. Results show that the isobaric interferences from Sb and Te are effectively eliminated in the LRIMS method.

A tin sample extracted from the spent fuel reprocessing solution is measured by the LRIMS method. The ratio of ^{126}Sn and $^{121\text{m}}\text{Sn}$ is measured as 49.78 without sample preparation chemistry procedure.

Conclusions

Based on the LRIMS device developed in the laboratory, a new attempt has been made for measurement of $^{126}\text{Sn}/^{121\text{m}}\text{Sn}$ ratio in the spent fuel reprocessing solution. The isobaric interferences are eliminated in this method. The results show that LRIMS may be the most appropriate technique to correctly determine isotope ratios if isobaric interferences exist from other elements in the sample.

Novel Aspect

Based on the LRIMS device developed in the laboratory, the auto-ionization energy spectrum of tin atom is obtained, and an efficient three-color-three-photon resonance ionization scheme of Sn is confirmed. A LRIMS method of measurement of $^{126}\text{Sn}/^{121\text{m}}\text{Sn}$ ratio in the sample extracted from the spent fuel reprocessing solution is proposed, in which the isobaric interferences caused by ^{126}Te and ^{121}Sb are eliminated.

References

1. Banghong Liang, Jinsong Zhang, Yicai Li et al., Journal of Chinese Mass Spectrometry Society, 33, 349-352 (2012).
2. Xiaopan Shen, Zhiming Li, Wenliang Wang, et al., Chinese Journal of Analytical Chemistry, 45, 342-347 (2017).
3. Brennecke G A, Borg L E, Romaniello S J et al., Geochimica et Cosmochimica Acta, 201, 331-344 (2017).
4. Oberli F, Gartenmann P, Meier M et al., International Journal of Mass Spectrometry., 184, 145-152 (1999).
5. Shen H, Jiang S, He M, et al., Nuclear Instruments & methods in Physics Research, 269, 392-395 (2011).
6. Fassett J D, Moore L J, Travis J C et al., Science, 230, 262-267(1985).

TP-98 / IMPROVING COVERAGE IN TRANSMEMBRANE DOMAINS OF INTEGRAL MEMBRANE PROTEINS WITH ULTRA-VIOLET PHOTODISSOCIATION

Julian Whitelegge

University of California Los Angeles, The Pasarow Mass Spectrometry Laboratory, Los Angeles

Keywords: membrane protein, FTMS, electron-capture dissociation, electron-transfer dissociation

Introduction

The most common structural feature of membrane proteins is the bilayer embedded alpha helix. Mass spectrometric coverage during dissociation experiments is frequently limited due to lack of charged

residues and retention of structure in the gas phase. In this presentation we explore new strategies for improving coverage in transmembrane domains using FTICR and other high-resolution instruments, including ultra-violet photodissociation (UVPD).

Methods

Samples were typically dissolved in formic acid immediately prior to separation by reverse-phase chromatography on polymeric stationary phase at 40°C (PLRP/S, Agilent) or size-exclusion in chloroform/methanol/aqueous formic acid (4/4/1). Samples were selected for high-resolution top-down mass spectrometry directed by low-resolution MS data and analyzed by static infusion freeing us from the constraints of chromatographic timescales.

Results

Samples were subjected to positive-ion nano-electrospray ionization mass spectrometry using either a 7 Tesla FTICR instrument (LTQFT Ultra with ECD and IRMPD) or an orbitrap (Lumos ETD with a solid state UV laser (213 nm) for UVPD). Addition of supercharging reagents shifted the predominant molecular ion to higher charge states, improving yield of product ions in ECD though it is unclear whether this was due to increased abundance of the precursor or denaturation by the supercharging reagent. We tested dissociation of the two transmembrane helix c-subunit of the ATP-synthase. The 7997.36 Da (monoisotopic) protein isolated from *S. oleracea* retained its initiating formyl-Met residue and yielded either 4- or 5-charge ions (m/z 2001 and 1334 respectively) both of which were analyzed by static nano-electrospray to yield a precursor intensity of $1e6$. These ions were subjected to various dissociation techniques including higher energy collisional (HCD), electron transfer (ETD), the binary mixture these two (ETHcD) for comparison with UVPD.

Conclusions

Modern dissociation strategies are giving excellent coverage in transmembrane domains of integral membrane proteins. While ETHcD enabled coverage above 50% for single charge states compared with CAD, HCD and ETD, UVPD gave 70-80% coverage. Summing UVPD data for two charge states gave close to 100% coverage.

Novel Aspect

UVPD and ETHcD appear especially effective for detailed structural characterization of the transmembrane domains of integral membrane proteins.

TP-99 / STUDY FOR MASS RESOLUTION OF QUADRUPOLE MASS SPECTROMETERS ACCORDING TO THE ROD ELECTRODE ALIGNMENT

Kiyomi Yoshinari⁽¹⁾ - ***Yasushi Terui***⁽²⁾

Hitachi, Ltd. Research & Development Group, Center for Technology Innovation - Controls, 1-1, Omika-cho 7 chome, Hitachi-shi, Ibaraki-ken (1) - Hitachi High-Technologies Corporation, Science & Medical Systems Business Group, 882, Ichige, Hitachinaka-shi, Ibaraki-ken (2)

Keywords: Quadrupole Mass Spectrometers, mass resolution, multipole electric field, rod electrode alignment, BEM

Introduction

Radio frequency (RF) quadrupole mass spectrometers (QMS) have come to be used in various fields because of their compactness. Ionized specimens are separated while passing through a QMS in accordance with their m/z values. The mass resolution of QMSs must be raised to identify species with high accuracy. The mass resolution is related to the alignment of a QMS's four rod electrodes [1].

Methods

We developed our original simulation program PISA-QMS, in which ion motions in an RF electric field can be calculated using the BEM [2,3] for an electric field and using the Shanks method for ion motion [4] within 0.01% error. PISA-QMS was applied to a QMS varying r/r_0 from 1.063 to 1.145, where r and r_0 indicate the radius of one rod and the inscribed circle between four rods. The relationship between mass resolution and r/r_0 values was investigated using PISA-QMS.

Results

The relationship between mass resolution (FWHM) and r/r_0 values showed a range of r/r_0 values with high mass resolution ($\text{FWHM} \leq 0.5$). The range of r/r_0 values giving high resolution ($\text{FWHM} \leq 0.5$) was $1.124 \leq r/r_0 \leq 1.13$. To clarify the mechanism of high mass resolution in such r/r_0 values, we calculated the electric fields between four rod electrodes and expanded them into multipole components ($n=1-6$). According to the values, the generation % of quadrupole components ($n=2$) and others ($n \neq 2$) changed. When r/r_0 was 1.127, the generation % ($n=2$) became almost 100% because the other multipole components totally canceled each other. The values of r/r_0 giving high resolution ($\text{FWHM} \leq 0.5$) show good agreement with the range of r/r_0 giving a high % of quadrupole component generation ($\geq 99.8\%$). High resolution could be obtained even when the rod electrode surfaces were not hyperboloid by generating a quadrupole electric field component greater than 99.8% and by canceling other components to no less than 0.2%.

Conclusions

An electric field, ion motions, and mass spectra were calculated by applying our simulation program PISA-QMS to QMS varying r/r_0 values. The results showed the mechanism of mass resolution is determined by quadrupole components % among the generated electric field. High resolution could be obtained even when the rod electrode surfaces were not hyperboloid by generating an electric field so as to make the quadrupole component $\geq 99.8\%$ and the other components $\leq 0.2\%$.

Novel Aspect

High resolution can be obtained by generating an electric field so as to make the quadrupole component $\geq 99.8\%$ and the other components $\leq 0.2\%$.

References

1. R. E. March and R. J. Hughes, Chemical Analysis Series Vol. 102, 42 (1989).
2. K. Yoshinari, et al., Nucl. Instrum. methods Phys. Res. Sect. A 363, 429(1995).
3. M. Koizumi and M. Utamura, Computational Mechanics 7, 183(1991).
4. E. B. Shanks, NASA TN D-2920 (1965).

TP-100 / SWEET HRMS WORKFLOW STRATEGY FOR SCREENING UNHEALTHY SUBSTANCES IN FOOD

Susanne Ekroth - Malin Johansson

National Food Agency, Chemistry Department, Uppsala, Sweden

Keywords: SweEt multiresidue method, Screening, Orbitrap mass analyzer, Compound Discoverer

Introduction

As more and more new food compounds are being developed and deployed around the world, contaminants not in target lists can appear in products through environmental or manufacturing contamination or misuse. This makes it essential that food safety laboratories stay vigilant for an ever-increasing number of analytes to screen for, and for new, unexpected and unwanted compounds[1] [2].

At the moment, the LC triple quad MS and GC triple quad MS systems are the workhorses of the pesticide laboratory. Most labs in Europe still use these technologies more than any other. But there's a trend toward implementing high-resolution accurate mass (HRAM) technology, because analysts are no longer satisfied with just finding the expected compounds in the samples; now, the aim is to identify everything (or at least as much as possible). In other words, there is a steady move away from targeted analysis. Screening allows us to discover what is in the samples with less bias – at odds with simply looking for what we already suspect is present.

This poster will show an efficient workflow for a non-targeted acquisition MS-mode for different compound groups and demonstrate a next-generation workflow for more efficient control of unexpected and unwanted compounds in food. The Swedish Ethyl acetate multiresidue method (SweEt), coupled with a high-resolution mass spectrometer offers:

The ability to screen for several compound groups at the same time, in the same sample

Advanced software tools for screening evaluation and identification

A validated and accredited method to screen for around 500 pesticides in fruits and vegetables

Methods

The main extraction techniques used today are making high-moisture food sample preparation faster and more straightforward. SweEt [3] (Swedish Ethyl acetate), is amenable to a wide range of matrices and offer a number of advantages over the traditional sample preparation techniques. With a range of polar and non-polar contaminants potentially present in food samples, it is essential that the sample preparation approaches are able to extract a very broad range of analytes. However, it can be challenging to develop methods that are both comprehensive as well as convenient to use. This presentation will describe the workflow for SweEt with around 550 pesticide analytes in targeted mode and the screening workflow for semi/non-targeted mode.

Results

With a generic method it is possible to screen and identify several compound groups in the same sample. The poster will present examples of findings of aflatoxin, benzophenone and pesticides.

Conclusions

Pesticides alone account for around 1000 analytes and can be found in different food commodities. Meanwhile other types of unwanted analytes can also be found such as mycotoxins, alkaloids, food additives, veterinary drugs and food contact materials. The broad variety of compound groups has limitations regarding extractability with different solvents. Therefore it is important to have a generic method and if possible to extract suspect samples in various ways. [4]

Novel Aspect

This poster presents an efficient screening workflow which can make laboratories better equipped to manage what the future may hold. There is a need to stay alert and anticipate what's beyond the conventional and expected findings, for example food frauds, contaminations or simply misuse of substances.

References

1. How to screen and identify unexpected and unwanted compounds in food. Thermo Fisher Scientific Application Note 65144 (2018).
2. How to screen and identify unexpected and unwanted compounds in food with SweEt multiresidue method. Webinar Thermo Fisher Scientific 19 Apr 2018.
3. Analysis in foods with ethyl acetate extraction using gas and liquid chromatography with tandem mass spectrometric determination (NMKL 195, 2013), NMKL/Publications/Methods/Methods-Chemical/NMKL 195 (<http://www.nmkl.org>)
4. Variety: the Spice of Life The Unofficial Guide to EPRW 2018, the Analytical Scientist https://theanalyticalscientist.com/fileadmin/tas/issues/18-05/0518_Thermo_EPRW_PR.pdf

TP-101 / SIMPLE WORKFLOW FOR PROTEOFORM SEPARATION AND ANALYSIS USING 2D GEL ELECTROPHORESIS-HIGH RESOLUTION MASS SPECTROMETRY

Manasi Gaikwad - Benjamin Dreyer - Hartmut Schlüter

University Medical Center Hamburg Eppendorf, Mass Spectrometric Proteomics Institute for Clinical Chemistry and Laboratory Medicine, Hamburg, Germany

Keywords: 2D gel electrophoresis, mass spectrometry, proteoforms, post translational modification

Introduction

In the early days of proteomics, it was assumed that each spot on a 2-dimensional electrophoresis (2DE) gel represent a single protein (1) due to the limitations of MALDI peptide mass fingerprinting for protein identification (2). The aim of this study was to answer the question how many proteins can be identified within one 2DE spot applying state of the art LC-MSMS techniques and if additional information about proteoforms can be obtained due to separation of proteins by 2DE.

Methods

Different 2DE-spots from a differential proteomics approach comparing human breast cancer cells before (control) and after perturbation (inhibition of mannosidase 1 with kifunensine) were analyzed. 60 spots were selected for further analysis after comparing the spot patterns of the 2DE gels. Selected spots were cut out and the proteins were proteolytically cleaved by trypsin. The peptides were analyzed by nano RPLC-ESI-MS/MS on a Q Exactive mass spectrometer.

Results

Thousands of proteins and proteoforms were separated on the 2DE gel seen visually on the gel as multiple spots. The 60 spots subjected to bottom up LC-MS/MS analysis resulted in having greater insight into the number of proteins behind each spot and the distribution of proteoforms across 2DE gels. On average 125 proteins were identified per spot of the control 2D gel. Using this strategy, proteins could be identified unambiguously. About 42 % of the identified proteins were found in more than one gel spot which suggest that the proteins have varied in the post translational modifications (PTMs). Putative annexin A-2 like protein, desmoplakin and junction plakoglobin were proteins which are most commonly found in multiple spots. These multiple protein spots represent proteoforms having varied isoelectric points and molecular masses but are derived from same parent gene.

Conclusions

The identification of proteins with state of the art LC-MS/MS from stained 2DE spots, which differ in their intensity in the comparison control vs. perturbed state gives information of proteoform changes associated with response to the perturbation. Combination of 2D gel separation and LC-MS/MS offers great potential in the proteome research focusing on proteoforms.

Novel Aspect

The workflow of 2DE followed by high resolution mass spectrometry represents a simple yet effective approach for identifying changes in the proteome on the level of proteoforms.

References

1. Person M., Shen J., Traner A., Hensley S., Lo H., Abbruzzese J. Li D., *Journal of BiomolecularTechniques*, 2006, 17, 2, 145–156.
2. Kleinert P., Kuster T., Arnold D., Jaeken J., Heizmann C., Troxler H.,*Proteomics*, 2007, 7, 15–22.

TP-102 / RAPID AND EFFECTIVE MULTI-CLASS ANALYSIS OF ANTIBIOTICS IN FEEDINGSTUFFS AT CARRY-OVER LEVEL BY LC-MS/MS TECHNIQUES

Ilaria Di Marco Pisciotto⁽¹⁾ - *Luigi Giannetti*⁽²⁾ - *Bruno Neri*⁽²⁾ - *Maurizio Fiori*⁽³⁾ - *Cinzia Civitareale*⁽⁴⁾ - *Paolo Stacchini*⁽⁴⁾ - *Pasquale Gallo*⁽¹⁾

Istituto Zooprofilattico Sperimentale del Mezzogiorno, Department of Chemistry, Portici (NA) (1) - Istituto Zooprofilattico Sperimentale del Lazio e della Toscana, Department of Chemistry, Rome (2) - Istituto Superiore di Sanità, National Reference Laboratory for additives in feed, Rome (3) - Istituto Superiore di Sanità, Department of Food Safety and Veterinary Public Health, Rome (4)

Keywords: antibiotic, mass spectrometry, carry-over, feed.

Introduction

The rough clean up of pipelines after medicated feeding stuffs production can result in a cross-contamination of feed, known as “unavoidable carry-over for non target species”. This can expose farmed animals to undesirable administration of antibiotics, increasing the risk of induced resistance to microbials. Therefore, LC-MS/MS based multi-class methods have been developed and compared to determine up to 38 antibiotics in feed for target and untargeted species.

Methods

Purification of up to 38 antibiotics from feed (tetracyclines, sulphonamides, quinolones, macrolides, pleuromutilins, diaminopyrimidines, streptogramins) was performed by solid-liquid extraction. Determination was carried out by three different analytical systems: LC/ESI-QTRAP-MS/MS, UHPLC-Orbitrap-MS and UPLC-Xevo-TQ-MS. In all cases, chromatographic separation was performed on reversed phase columns, using linear gradient elution.

Results

The three test methods were in-house validated according to requirements of Commission Regulation 882/2004/EC [1] and Decision 657/2002/EC [2], evaluating the analytical performance parameters specificity, linearity, limit of quantification (LOQ), limit of detection (LOD), trueness, precision and ruggedness for each analyte. Specificity was assessed analyzing uncontaminated feed samples for poultry, cattle, pig, rabbit, sheep, goat and equine; no matrix interference was observed in the chromatograms at analyte retention times. Trueness as mean recoveries (%), and precision as RSD (%), were evaluated for each method by analyzing blank feed samples spiked at least at three concentration levels and were, respectively, in the range 76.3-107.7% and 1.6-22.8%. Linearity of each detector was verified by the correlation coefficient ($r^2 \geq 0.950$) of linear regression calibration curves of each analyte. The LOQ for all the antibiotics ranged from 0.010 to 0.100 mg/kg depending on the equipment; likely, LOD ranged from 0.002 to 0.030 mg/kg.

Conclusions

The multi-class and multi-drug methods we developed are all based on a rapid and simple cleanup, an effective chromatographic separation and a highly sensitive detection, allowing for simultaneous quantification and unambiguous identification of up to 38 antibiotics down to 0.010 mg/kg. These analytical methods fit with the requirements of the European Committee for the official control of antibiotics in feed at carry-over levels.

Novel Aspect

These analytical methods can replace single-class methods resulting in an optimization of work, costs and staff organization and in a more effective control activity.

References

1. The European Commission, Official Journal of European Union, L 165, 1-141 (2004).
2. The European Commission, Official Journal of European Union, L 221, 8-36 (2002).

TP-103 / DETERMINATION OF ORGANIC ARSENIC COMPOUNDS BY MEANS OF HIGH RESOLUTION MOLECULAR MS/MS

Ronald Glabonjat - Kenneth Jensen - Nikolaus Guttenberger - Kevin Francesconi

University of Graz, Analytical Chemistry, Graz, Austria

Keywords: Organic Arsenic Species, High Resolution ES-MS/MS, Metabolomics, Fragmentation, Orbitrap.

Introduction

Arsenic is unique among trace elements by exhibiting a rich organic chemistry spanning over 100 naturally occurring organoarsenic species in the environment. Organoarsenic species range from simple methylated arsenic acids through to arsenic-bound sugars and even more complex arsenic-containing lipids [1]. Here we report our most recent results employing high resolution mass spectrometry (HR-MS) to identify novel environmental arsenic compounds.

Methods

We characterize the complex organoarsenical compositions in various environmental samples by means of positive ionisation HR electrospray Orbitrap MS/MS determination after separation of individual species by reversed-phase HPLC. The importance of tandem-MS as a powerful tool for the arsenic specific identification of unknown compounds is demonstrated and factors influencing the selectivity and sensitivity of the method discussed in detail.

Results

The composition of the naturally occurring organoarsenic compounds are revealed by molecular HR-MS/MS measurements, whereby most of them show characteristic fragmentation. Among these, the most significant fragments $[M+H]^+$ have m/z of 104.9680 ($C_2H_6As^+$) and 122.9786 ($C_2H_8AsO^+$) originating from a common polar dimethylarsinoyl head group $(CH_3)_2OAs$. Furthermore, there is a variety of other fragmentation products supporting the identification of structural features such as arseno-ribosides with m/z 237.0102 ($C_7H_{14}AsO_4^+$), or arseno-phospholipids m/z 483.0607 ($C_{13}H_{29}PAsO_{12}^+$). In dexterous combination with non-arsenic containing fragments such as *e.g.* 184.0733 ($C_5H_{15}O_4PN^+$) or corresponding neutral losses *e.g.* m/z 183.0660, we also show simple ways to identify even more complex arsenic compounds such as cholines or ethanol-amines.

Conclusions

In our work we provide a sophisticated but still simple way for the determination of a wide range of organoarsenicals using high resolution mass spectrometry. Our method is not dependent on standard-compounds, and therefore has potential to open up the field for other scientists and assist in further research in this rapidly expanding field of environmental research.

Novel Aspect

We present an approach to species independent, arsenic-specific determination without the necessity of an element-specific detector by using HR-ESMS/MS and monitor characteristic organoarsenic fragments.

References

1. Glabonjat R. A., Raber G., Jensen K. B., Guttenberger N., Zangger K., Francesconi K. A., *Angew. Chem. Int. Ed.*, 56, 11963–11965 (2017).

TP-104 / LINEAR DYNAMIC RANGE IMPROVEMENT IN TOF-MS/MS MODE ON X500 QTOF SYSTEM

Feng Zhong - Suyu Liu - Wen Jin - Doug Simmons - Nic Bloomfield - Robert Haufler

Sciex, Research, Concord, Canada

Keywords: HRMS, TOFMS, LDR, QTOF, Quantitation

Introduction

Linear Dynamic Range for MS-MS mode was improved by employing beam modulation. This technology was previously used in MS mode, but has now also been developed for MS-MS mode. This allowed the extension of LDR by increasing the upper limit of quantitation, with no changes to the lower limit of quantitation.

Methods

17 standards were prepared by making 3-fold serial dilutions of 100 µg/mL Alprazolam/Atorvastatin in 50/50 acetonitrile/water containing 10 ng/mL Alprazolam-D₅/Atorvastatin-D₅, respectively. The analyte concentration covered the range of 0.00077400 to 33,300 ng/mL.

Experiments ran alternating between modulation on and modulation off acquiring spectra in TOF-MS/MS mode.

Results

In each case, the LLOQ was not affected by EDR, and the improvement of LDR is due to extension of ULOQ to higher concentration samples as expected. Accumulation times down to 25 ms were tested and did not show a significant impact on LDR. The average increase in the LDR was about 100 fold, from about 2.5 orders up to about 4-4.5 orders.

The EDR feature for MS/MS improves the LDR performance for tested compounds and workflows by extending the ULOQ to higher concentration without effects on LLOQ.

Shorter accumulation, 25 ms and 35 ms in MS/MS mode has not caused any loss of LDR performance comparing with 50 ms accumulation time.

References

1. Guidance for Industry Bioanalytical Method Validation, FDA. May 2018.

TP-105 / APPLYING A NOVEL COMPONENT DETECTION ALGORITHM TO HELP ACCELERATE METABOLOMICS DISCOVERY WORKFLOWS

Christine Hinz⁽¹⁾ - ***Emily G Armitage***⁽²⁾ - ***Kirsten Hobby***⁽²⁾ - ***Neil Loftus***⁽²⁾ - ***Nicola Gray***⁽³⁾ - ***Manuel Y Schar***⁽³⁾
- ***Jeremy P E Spencer***⁽³⁾

Shimadzu UK Limited, LCMS technical department, Milton Keynes, United Kingdom (1) - Shimadzu Corporation, Research and Development, Manchester, United Kingdom (2) - University of Reading, Department of Food and Nutritional Sciences, School of Chemistry, Food and Pharmacy, Reading, United Kingdom (3)

Keywords: Component detection, Peak picking, Metabolomics, Untargeted analysis, Accurate mass

Introduction

In untargeted metabolomics workflows, data mining techniques are designed to identify differences between metabolite profiles from control and test groups, relating to specific clinical challenges or states. To accelerate metabolite detection we've developed the novel algorithm 'FIND' to locate components using high mass accuracy LC-MS TOF data. Here, it is applied to assess the impact of phenolic acids and avenanthramides in the diets of adult human volunteers.

Methods

Blood was provided before and after a 4-week intervention consuming a daily moderate dose of phenolic acids and avenanthramides in oats or a phenolic-low control. Extracted plasma was analysed using high mass accuracy LC-MS/MS TOF (Shimadzu Corporation, Japan). The FIND algorithm was applied to raw data

with only an estimate of typical chromatographic peak width (FWHM) used. A collection of inbuilt logic is used to correctly identify masses of the unknown components.

Results

Untargeted discovery-based metabolomics workflows are designed to detect components from complex sample matrices. To help accelerate component detection in metabolomics workflows a novel algorithm has been designed to locate ion signals that are covariant and related in context (charge-state, isotope distribution, adducts, multimers) so that only a true monoisotopic mass of the component is reported.

Using a food-omics study as a model system, the component detection algorithm was applied to analysis of metabolite profiles from a human volunteer study exploring the impact of phenolic acids and avenanthramides in oats. In this study, the component detection algorithm was used to locate individual components in raw data files and to help verify metabolite detection between each test group following dietary exposure. The application of the component detection algorithm is simple, takes seconds to process, and does not require a complex method with many user parameters.

Conclusions

Adult volunteers provided blood before and after a 4-week intervention consuming a daily moderate dose of phenolic acids and avenanthramides in oats or a phenolic-low control diet in a crossover design. LC-MS/MS TOF was used to analyse extracted plasma and the FIND algorithm was applied to raw data to correctly identify the masses of the unknown components. A collection of inbuilt logic correctly identifies monoisotopic masses and related ions are grouped.

Novel Aspect

Applying a novel algorithm to high mass accuracy LC-MS TOF data to locate individual components in complex samples.

TP-106 / COMPREHENSIVE ANALYSIS OF FLAME RETARDANTS AND OTHER HALOGENATED CONTAMINANTS IN DUST FROM AN E-WASTE RECYCLING CENTRE

Alin Ionas⁽¹⁾ - *Karl Jobst*⁽²⁾ - *Lisa Melymuk*⁽¹⁾ - *David Bowman*⁽³⁾ - *Miriam Diamond*⁽⁴⁾ - *Victoria Arrandale*⁽⁵⁾ - *Marta Venier*⁽⁶⁾ - *Liisa Jantunen*⁽⁷⁾ - *Eric Reiner*⁽³⁾ - *Jana Klánová*⁽¹⁾

Masaryk University, RECETOX, Brno, Czech Republic (1) - McMaster University, Chemistry department, Hamilton, Canada (2) - Ontario Ministry of the Environment, Laboratory Services Branch, Toronto, Canada (3) - University of Toronto, Department of Earth Sciences, Toronto, Canada (4) - University of Toronto, Dalla Lana School of Public Health, Toronto, Canada (5) - Indiana University Bloomington, School of Public and Environmental Affairs, Bloomington, United States (6) - Environment and Climate Change Canada, Centre for Atmospheric Research Experiments, Egbert, Canada (7)

Keywords: halogenated chemicals, GC × GC-HRMS, compound identification, GC-QTOFMS, dust analysis

Introduction

Consumer electronics typically contain halogenated flame retardant chemicals to comply with flammability standards. However, various negative health effects have been associated with the aforementioned compounds [1].

The focus of this study was analysis of selected dust samples from an e-waste dismantling facility in Toronto, Canada. The objective was to determine any halogenated additives present in the e-waste, and especially additives which may not yet be described in the literature.

Methods

The samples were extracted using solvent-liquid extraction and ultrasound-assisted extraction and did not undergo destructive clean-up.

The extracts were injected on a Waters Xevo G2-XS quadrupole time-of-flight mass spectrometer coupled to an Agilent 7890B GC with atmospheric pressure chemical ionization (APCI) in two configurations:

- “Fast GC”, using an Agilent DB5-HT 15 m × 0.25 mm internal diameter × 0.1 μm film thickness; a (He) carrier gas flow of 3 mL/min and a 10 min run time.
- Comprehensive two-dimensional gas chromatography (GC × GC), using a method with a slow temperature gradient and an Agilent DB-5 60 m × 0.25 mm × 0.25 μm as the first dimension column and a Restek Rtx-17SIL (1 m × 0.15 mm × 0.15 μm) second dimension column. The 2° oven was offset +15°C from the 1° oven.

Results

In the fast GC screening, the halogenated compounds tentatively detected were mainly flame retardants, such as polybrominated diphenyl ethers (and especially BDE-209, which was highly concentrated in some samples), TBBPA and BEHTBP (both present in very high amounts) and EHTBB. Polychlorinated biphenyls were also detected in the samples, but in lower amounts than the flame retardants.

The GC×GC data was further investigated for the presence of target compounds, by extraction of a selected ion chromatogram (SIC) for all target masses. Non-target compounds were also identified on the basis of visual inspection of Kendrick Mass Defect (KMD) plots, constructed using-H/+Cl substitution (34 Da)[2].

The list of halogenated isotopic clusters generated from the 1D-GC data was searched in the GC×GC data (GC Image HRMS v2.5) where, due to the higher separation of peaks, co-elutions with other analytes are less likely. Consequently, the spectra of the target and non-target compounds of interest were “pure”. Structure assignment was aided by NIST library searches. The (soft) nature of the ionization, *i.e.* APCI in the positive mode, typically provides the molecular ion peak as the base peak with few fragments. Despite this, the matches obtained with the NIST library (containing electron ionisation spectra) were quite good, because the search algorithm places more weight on the matches of higher molecular mass peaks/fragments.

Using this technique, a methoxy derivative of TBBPA was tentatively identified, as were a few dibromophenyl and tribromophenyl isomers. A few non-halogenated compounds were detected as well, mainly organophosphate esters. Among these, the tricresyl phosphates and other triaryl phosphates were in particularly high amounts.

Conclusions

The dust from an e-waste dismantling site contained a wide variety of flame retardants, so dust collected from these types of sites is a good indicator of which flame retardants are contained in various types of consumer electronics.

The approach presented in this study has the benefits of providing quick information through the Fast GC screening. The enhanced separation of GC×GC provided in-depth information (*viz.* isobarically pure mass spectra and multidimensional GC retention behavior) at the cost of analysis time.

Novel Aspect

To the authors’ knowledge, this is the first study using both a Fast GC screening and GC × GC for identification and confirmation of chemicals identified in e-waste dust samples.

References

1. Lyche J. L., Rosseland C., Berge G., Polder A. *Environ. Int.* 74, 170–180, (2015).
2. Taguchi V. Y., Nieckarz R. J., Clement R. E., Krolik S., Williams R. J. *Am. Soc. Mass Spectrom.* 21 (11), 1918–1921, (2010).

TP-107 / METABOLOMIC AND N-GLYCAN CHARACTERIZATION OF SKIN FIBROBLASTS FROM AGED AND YOUNG DONORS

***Anna Karen Carrasco Laserna*⁽¹⁾ - *Bao Hui Ng*⁽¹⁾ - *Jeyakumar Masilamani*⁽²⁾ - *Toan Thang Phan*⁽²⁾ - *Fong Yau Sam Li*⁽¹⁾**

National University of Singapore, Department of Chemistry, Singapore (1) - CellResearch Corporation, National University of Singapore, Department of Surgery, Singapore (2)

Keywords: skin fibroblasts, aging, metabolomics, *N*-glycans

Introduction

Fibroblasts are cells that play an important role in wound healing. They are responsible for excreting extracellular matrix at the site of injury and in their myofibroblast form, cause wound contraction and deposit collagen[1]. Delayed wound healing is commonly observed with aging and a major contributing factor is the dysfunction of the fibroblasts [2]. In this study, we explore the metabolic and protein glycosylation changes in skin fibroblasts with aging.

Methods

Skin fibroblasts were collected from young (YSF) and aged (ASF) donors, cultured to $\sim 1 \times 10^6$ cells, and extracted by sonication with 80% methanol. The supernatant was obtained and analyzed using reverse phase (RP) and HILIC LC-MS in positive and negative ESI modes. Precipitated proteins were separated from cell debris, deglycosylation was performed and the *N*-glycans were labeled using 2-AB. The samples were also analyzed using RP and HILIC LC-MS in positive ESI mode.

Results

Initial multivariate analysis using Principal Components Analysis (PCA) of metabolomic LC-MS data from positive ESI HILIC and RP modes showed distinct differentiation of ASF and YSF. Orthogonal Projection to Latent Squares- Discriminant Analysis (OPLS-DA) was also performed to identify features that are significantly contributing to the differentiation of the ASF and YSF. Identification of differentiating features revealed that certain phospholipids were lower in ASF samples while acylcarnitines were found to be in higher levels in ASF. This could be indicative of higher rate of lipolysis and reduced mitochondrial fatty acid oxidation in ASF. *N*-glycan LC-MS analysis also revealed differences in the glycan profile for the ASF and YSF samples. Identification of the significantly differentiating glycans is on-going.

Conclusions

Preliminary results revealed significant differences in the metabolic and *N*-glycan profile of skin fibroblasts from young and aged donors. Skin fibroblasts from aged donors were found to have lower phospholipid levels and higher acylcarnitine levels indicating higher rate of lipolysis and reduced mitochondrial fatty acid oxidation function.

Novel Aspect

We characterized the metabolic and *N*-glycan profiles of skin fibroblasts from a single set of sample extraction revealing significant changes in skin fibroblasts with aging.

References

1. Werner S., Grose R., *Physiol Rev*, 83, 835-870 (2003).
2. Sgonc R., Gruber J., *Gerontology*, 59, 159-164 (2013).

TP-108 / LEAN APPROACH OF AUTOMATED SAMPLE PREPARATION FOR LC-MS OF INTACT PROTEINS ON A LIQUID HANDLING ROBOT

Katharina List - Dorothee Wetzel - Martin Pauers

Boehringer Ingelheim GmbH & Co. KG, Biologicals Development, Biberach, Germany

Keywords: Automation, Sample Preparation, Intact Protein, Handling Robot, Biotherapeutics

Introduction

The preparation of samples for LC-MS analyses is elaborate and time consuming. Due to more routine utilization of LC-MS in biopharmaceuticals development, the amount of samples rises, raising the need to automate sample preparation. An automated process will enable faster sample preparation and higher reproducibility and repeatability of the sample.

Methods

A manual sample preparation protocol, consisting of dilution steps, was transferred to a robotic liquid handling system (Tecan, Fluent). An internal standard was used and added into the samples before the dilution steps for protein concentration determination by LC-MS (Waters Acquity UPLC coupled to a Waters SYNAPT G2-S). Additionally the protein concentration was determined by UV-photometer.

Results

LC-MS data confirms the UV data and shows that automated sample preparation presents higher quality in repeatability, reproducibility and comparability in contrast to manual sample preparation. Regardless of the initial concentration of the samples the robotic system is able to generate dilutions fitting more precisely to the target concentration. Furthermore it manages to reproduce the same concentration of a sample that was diluted earlier with only minimal deviation and still more exactly as if diluted manually. In addition the concentrations of the automated prepared dilutions are more accurate from sample to sample meaning an increase of comparability between the samples compared to manual sample preparation. The robotic liquid handling system from Tecan is capable of handling complex and high numbers of protein samples.

Conclusions

The sample preparation protocol was successfully transferred to a robotic liquid handling system. LC-MS data demonstrated that automated sample preparation is highly reproducible, repeatable and comparable. Further complex sample preparation, such as peptide mapping protocol, is planned to be transferred to the robot.

Novel Aspect

Reproducible implementation of automated sample preparation on a robotic liquid handling system was achieved.

TP-109 / COMBINATION OF CHEMICAL IONIZATION (CI) AND LOW-ENERGY ELECTRON IONIZATION WITH HIGH-RESOLUTION GC/Q-TOFMS

Viorica Lopez-Avila⁽¹⁾ - *Sofia Nieto*⁽¹⁾ - *Harry Prest*⁽¹⁾ - *Jeffrey Kernan*⁽¹⁾ - *George Yefchak*⁽¹⁾ - *Robert Clark*⁽¹⁾ - *Nathan Eno*⁽¹⁾ - *Jim Oppenheimer*⁽¹⁾ - *Bill Russ*⁽¹⁾ - *Nicola Cimino*⁽²⁾

Agilent Technologies, Mass Spectrometry, Santa Clara, United States (1) - Agilent Technologies, EMEA, Rome, Italy (2)

Keywords: chemical ionization, low-energy electron ionization, gas chromatography, Q-TOF mass spectrometry

Introduction

GC/Q-TOF mass spectrometry (MS) with low-energy electron ionization (EI) and chemical ionization (CI) have developed into versatile analytical tools for the compound identification due to the enhancement in the relative abundance of the molecular ions or their adducts, respectively, as compared to the standard EI. The latter technology offers additional advantages as the ionization process can be controlled with various reagent gasses. The use of combination of CI and GC/Q-TOF MS will be shown here with several compound classes including organochlorine and organophosphorus pesticides.

Methods

Experiments were done with an Agilent 7250 GC/Q-TOF equipped with a low-electron energy EI and an interchangeable CI source. The test compounds included: organochlorine (OCP) and organophosphorus (OPP) pesticides, chlorinated acids and a few conazoles. The GC separation was on a 30 m x 0.25 mm id x 0.25 μ m HP-5MS capillary column using He at 1.2 mL/min. The injector temperature and the MS interface were set at 280°C. Methane (99.99%) was used as reagent gas. For positive CI, the methane flow was set at 20% and for negative CI at 40%. For NCI, the source and the quadrupole temperatures were set at 150°C. For PCI, the source temperature was set at 280°C and the quadrupole temperature at 150°C. The spectral data were acquired at 5 Hz and the mass range was 50-650 m/z . 2H-Perfluoro-5,8-dimethyl-3,6,9-trioxadodecane (PFDTD) was used to tune the mass spectrometer in the CI mode. Data analysis was performed using Agilent MassHunter (MH) Qualitative Analysis software version B.08 as well as MH Quantitative Analysis software version B.09.

Results

Low-energy EI provided an increase in relative abundance of molecular ions for halogenated dienes like hexachlorobutadiene and hexachlorocyclopentadiene, halogenated aliphatic ethers, and some phthalate esters. The interchangeable prototype CI source was evaluated with selected OCPs and OPPs in negative mode using methane as reagent gas. Significant decrease in fragmentation was observed for the OCPs, OPPs, and the chlorinated acids as methyl esters, when using negative CI, although the molecular ion was not always the base peak for these compounds. Nonetheless, fragmentation in CI leads to larger m/z fragment ions that helps in both compound identification as well as targeted analysis. The NCI spectra of a series of 14 conazoles including aza-, eta-, propi- and difenoconazole will be presented. NCI has a different fragmentation mechanism than EI and with high resolution Q-TOF such differences between the two ionization mechanisms can be easily explained.

Conclusions

Benefits of the 7250 GC/Q-TOF system equipped with a low energy-capable EI source as well as an interchangeable prototype CI source were explored for targeted and untargeted analysis applications. Chemical ionization alone or in combination with low-energy EI and a high-resolution GC/Q-TOF provides new opportunities in compound identification.

Novel Aspect

Use of high-resolution GC/Q-TOF MS with interchangeable low-energy EI or prototype CI sources in positive or negative mode

TP-110 / LIQUID CHROMATOGRAPHY COUPLED TO QUADRUPOLE-ORBITRAP HIGH RESOLUTION MASS SPECTROMETRY BASED METHOD FOR TARGET ANALYSIS AND SUSPECT SCREENING OF NON-IONIC SURFACTANTS IN TEXTILES

Xin Luo - Zengyuan Niu - Xiwen Ye - Zhixu Tang - Li Zhang

Qingdao Customs, Inspection and Quarantine Technical Center, Qingdao, China

Keywords: Alkylphenolpolyethoxylates (APEOs), Alcohol polyethoxylates (AEOs), Orbitrap, High resolution mass spectrometry, Textiles

Introduction

Non-ionic surfactants such as alkylphenolpolyethoxylates(APEOs) and alcohol polyethoxylates (AEOs) are widely applied in textile industries. They are potential sources of environmental pollution. In recent years, APEOs and AEOs have been detected and reported frequently in environmental samples, such as water and sediments. In this study, we describe a high-throughput and sensitive method for textiles analysis, using liquid chromatography coupled to quadrupole-Orbitrap high resolution mass spectrometry(LC-Q-Orbitrap

HRMS), for the simultaneously quantitative analysis of 40 target APEO oligomers with reference standards and screening of 160 AEO oligomers without standards in textiles.

Methods

The APEOs contain nonylphenol ethoxylates (NPEOs) and octylphenol ethoxylates (OPEOs) with an EO number of ethylene oxide of 1-20, while AEOs focus on C11EOs-C18EOs with an EO number of 1-20. After ultrasonic extraction in methanol, the extract was directly separated using a core-shell CORTECS C18+ column and analyzed by Full MS/dd-MS2 (data dependent acquisition) scan in ESI positive mode. Two best sensitivity experimental conditions for APEOs with short EO chains (AP(EO)1-2) and long EO chains (AP(EO)3-20) were investigated, respectively.

Results

Most APEO oligomers had wide concentration ranges and the correlation coefficients (R^2) were higher than 0.999. The limit of quantitation (LOQ) values for NP(EO)3-20 oligomers ranges from 16.00 to 52.80 $\mu\text{g}/\text{kg}$ and for OP(EO)3-20 oligomers is from 2.40 to 8.00 $\mu\text{g}/\text{kg}$. LOQ for NP(EO)1 and NP(EO)2, OP(EO)1 and OP(EO)2 was 2.40 mg/kg and 0.24 mg/kg, 1.20 mg/kg and 0.16 mg/kg, respectively. The average recovery for each APEO oligomer in cotton and polyester matrix was between 78% and 110% at three spiked levels and the relative standard deviation (RSD %) was below 10%. Theoretical m/z and assignments for MS/MS product ions of AP(EO)1-20 were predicted, the fragmentation regularity was divided into 7 groups according to different structures and three typical ions were listed. At the same time, based on the similarity of APEOs and AEOs in structure as well as previous reports on these two kinds of non-ionic surfactants, three typical product ions of APEOs were selected for confirmation of AEOs without reference standards. Then for AEOs suspects, a HRMS compound database containing 160 AEO oligomers was built and several parameters such as exact m/z , isotopic patterns, predicted product ions and predicted retention time were used for screening and confirmation. The established method was successfully applied for analysis of 40 commercial textile samples. Compared with OPEOs, NPEOs, especially NP(EO)3-15 oligomers, were widely detected in samples and the total concentration ranged from 1.56 to 1376.31 mg/kg. AEOs were also found in most samples, among which C12-14, C16 and C18 compounds appeared more frequently and the EO chains mainly ranged from 3 to 15.

Conclusions

A high-throughput and sensitive LC-Q-Orbitrap HRMS method has been established for simultaneously quantitative analysis of 40 APEOs and screening of 160 AEOs in textiles. Under two optimized chromatographic and mass spectrometric conditions, all oligomers with EO chains from 1 to 20 can be analyzed with high sensitivity. The method has been successfully applied to analyze 40 textile samples. NP(EO)3-15 and AEO homologs of C12-14, C16 and C18 with EO chains ranging from 3 to 15 were frequently detected in samples. This method is rapid, high throughput and sensitive. Moreover, any other potential APEOs and AEOs oligomers interested can be added into the compound database and retrospective analysis of them is feasible for present method.

Novel Aspect

40 APEOs targets and 160 AEOs suspects in textiles were simultaneously analyzed. APEOs with both short and long EO chains can achieve high sensitivities. Exact m/z and fragmentation pathways of main APEOs fragments were predicted. Fragmentation regularities of APEOs were used for confirmation of AEOs suspects. Method was applied on 40 textile samples and oligomers detected were concluded.

TP-111 / COMPREHENSIVE ANALYSIS OF MULTI-CLASS HAZARDOUS CHEMICALS IN PLASTIC PRODUCTS USING HIGH RESOLUTION MASS SPECTROMETRY

Zengyuan Niu - Xin Luo - Xiwen Ye - Zhixu Tang - Li Zhang

Qingdao Customs, Inspection and Quarantine Technical Center, Qingdao, China

Keywords: multi-class, hazardous chemicals, plastic products, high resolution mass spectrometry

Introduction

With the increasing awareness of the safety of consumer goods, the relevant restrictions on the number of hazardous chemicals are increasing. Simultaneous detection for multi-class targets and non-targets is becoming a more popular and difficult topic.

Methods

Present work established a novel method for the comprehensive analysis of 41 chemical hazards (20 targets and 21 non-targets) in plastic products (including plastic toys, electrical and electronic plastic parts, polymeric food contact materials, etc.). Compounds studied here belong to four classes, mainly including phthalic acid esters (PAEs), organotins (OTs), perfluorochemicals (PFCs) and flame retardants (FRs), et al. According to the sample material, appropriate organic solvents were selected for ultrasonic extraction, and methanol was used to precipitate the polymer. After standing for some time, the extract was separated by the core-shell particle column and analyzed by the Orbitrap high resolution mass spectrometry.

Results

The results were confirmed by multidimensional identification of the parameters such as the m/z of the parent ion, the isotope ratio, and the exact m/z of fragments. Method validation for 20 compounds was evaluated. The limits of quantitation ranged from 2 mg/kg ~ 200 mg/kg, the recovery was 74% ~ 120% and the precision was less than 13%. Present method was applied for simultaneous detection both targets and non-targets in 34 plastic samples, among which di(2-ethylhexyl) phthalate (DEHP), dibutyl phthalate (DBP) and diisobutyl Phthalate (DIBP) were detected quantitatively and tetrabromobisphenol A (TBBPA) was detected qualitatively.

TP-112 / HIGH RESOLUTION MASS SPECTROMETRY BASED METHOD FOR SCREENING AND IDENTIFICATION OF 48 PROHIBITED DYES IN TEXTILES USING DATA-DEPENDENT ACQUISITION MODE

Zengyuan Niu - Xin Luo - Zhixu Tang - Xiwen Ye - Lianyi Qu

Qingdao Customs, Inspection and Quarantine Technical Center, Qingdao, China

Keywords: high resolution mass spectrometry; screening and identification; dyes; textiles; data-dependent acquisition mode

Introduction

Various prohibited dyes are often misused or even deliberately used for pursuit of economic interest in textile industries. Owing to their widespread use, toxicity on human health, as well as destruction to the ecological environment, some international authoritative consumer laws have strictly limited their use. Dyes used in textiles often contain synthetic intermediates, dispersants and other impurities, making their composition very complex. While traditional analysis usually focus on single group of compounds, and false positive identifications might be possible due to the low selectivity of methods.

Methods

The currently presented work firstly developed a high-throughput screening and identification method for simultaneous analysis of 48 prohibited dyes, belonging to 5 different groups, in textiles. After extraction in pyridine/water (1/1, V/V) water bath at 95°C, the extract was filtered with 0.22 μm membrane, then directly separated using a C18 column and analysed by Full MS/dd-MS2 (data-dependent acquisition) scan using electrospray ionization probe (ESI) operating in positive and negative mode. The retention time, exact

mass of each precursor and MS2 fragment, as well as the isotopic pattern of precursors were used for the simultaneous analysis of targeted dyes. Moreover, for the identification of complex components within certain dyes such as disperse blue 35, software of ACD/Labs were also used for the prediction of MS2 fragment assignments and calculation of the theoretical accurate mass, which greatly improve the identification reliability and accuracy of the presented method.

Results

All 48 dyes had wide concentration ranges and the correlation coefficients (R²) were higher than 0.99. The limit of quantitation (LOQ) values ranged from 0.01 mg/kg to 10 mg/kg. The average recovery for each analyte in representative cotton and polyester matrix was between 62% and 120% at three spiked levels and the relative standard deviation (RSD%) was below 15%. The established method was successfully applied for analysis of 40 commercial textile samples with various fiber types, among which disperse orange 37 was detected. This method is rapid, sensitive and accurate, which could be applied on the high throughput analysis of prohibited dyes in textiles.

TP-113 / QUANTITATION OF SERUM CREATININE IN MICE USING HYDROPHILIC INTERACTION LIQUID CHROMATOGRAPHY (HILIC) COUPLED WITH HIGH-RESOLUTION ACCURATE-MASS (HRMS) ORBITRAP MASS SPECTROMETER

Hidenori Takada⁽¹⁾ - *Yuko Ube*⁽²⁾ - *Mariko Sakai*⁽²⁾ - *Yoko Hasegawa*⁽¹⁾ - *Tetsuya Kitayama*⁽²⁾ - *Yuji Yamazaki*⁽³⁾ - *Tatsuo Uchimura*⁽¹⁾

Kyowa Hakko Kirin Co., Ltd, Translational Research Unit, R&D Division, Sunto-gun, Japan (1) - Kyowa Hakko Kirin Co., Ltd, Nephrology R&D Unit, R&D Division, Sunto-gun, Japan (2) - Kyowa Hakko Kirin Co., Ltd, Nephrology R&D Unit, R&D Division, Machida, Japan (3)

Keywords: Creatinine, HILIC, HRMS, Biomarker, Cyanoacetic acid

Introduction

Serum creatinine level (SCR) is widely used as a biomarker of renal function. The SCRs are lower in normal mice than in human, dogs and rats. Colorimetric methods (Jaffe method and enzymatic method) commonly used in laboratory and clinic are insufficient from perspectives of accuracy and sensitivity for determination of the mouse SCR. A highly accurate and sensitive analytical method using LC/MS/MS for determination of mouse SCR has been developed.

Methods

Creatinine was separated in HILIC mode using cyanoacetic acid as an acidic additive in LC mobile phase and detected by parallel reaction monitoring (PRM) in positive ESI mode using hybrid quadrupole-orbitrap mass spectrometer (Q Exactive Focus, Thermo Fisher). The SCR in mice with renal ischemia-reperfusion (I/R) and sham operation were determined using the novel HILIC/HRMS method and the colorimetric methods and compared with blood urea nitrogen (BUN).

Results

In the HILIC/HRMS method, cyanoacetic acid, whose decarboxylation occurs by heating at temperature higher than 160°C and produces acetonitrile with carbon dioxide, resolved the matrix effect with decrease in MS response of creatinine in mouse serum observed using formic acid as an acidic additive. The selectivity of creatine in the blank serum was evaluated using peak area ratio of quantifying ion mass (m/z 86.07) to qualifying ion mass (m/z 72.04) (Quant/Qual ratio), and the issue was minimal. The intra- and inter-assay accuracy and precision, and the results in stability testing (storage at -20°C for 95 days and freezing and thawing) were preferable. In I/R model and sham-operated mice, the correlation of SCR by the HILIC/HRMS method with BUN was better than those of SCRs by colorimetric methods (Jaffe method and

enzymatic method) with BUN. In particular, the deviations from the correlation of SCR by the HILIC/HRMS method were smaller at lower levels in normal mice or mice with mild renal failure than that of SCR by enzymatic method.

Conclusions

The novel analytical method for determination of serum creatinine in mice has been developed and utilizes separation in HILIC mode using cyanoacetic acid as an acidic additive of LC mobile phase and HRMS detection in positive ESI mode. Cyanoacetic acid resolved the matrix effect in MS response of creatinine in mouse serum. The developed method showed potentials of highly accurate determination for lower SCR (≤ 0.3 mg/dL) in normal mice or mice with mild renal failure.

Novel Aspect

The novel HILIC/HRMS method showed accurate and precise determination of mouse serum creatinine. Cyanoacetic acid showed a novel potential for LC/MS detection in ESI mode.

TP-114 / ULTRASENSITIVE DETECTION OF PROTEIN BIOMARKERS BY SIGNAL AMPLIFICATION MASS SPECTROMETRY

Yuning Wang - Ruijun Du - Liang Qiao - Baohong Liu

Department of Chemistry, Fudan University, Shanghai, China

Keywords: biomarker, membrane protein, mass spectrometry, signal amplification

Introduction

Proteins are essential biomolecule in revealing biological activities, while abnormal expression of proteins is often related to various diseases. [1, 2] However, protein biomarkers expressed in cells can be exceedingly rare. Moreover, unitary biomarker is usually expressed by more than one type of disease. [3, 4] Thus, sensitive and accurate detection of multiple biomarkers is indispensable for clinical diseases diagnosis.

Methods

Herein, we propose a novel mass spectrometry (MS) signal amplification method to detect low-abundance protein biomarkers using mass labels (MLs)-modified Au nanoparticles (AuNPs), based on our previous work. [5] AuNPs modified with specific antibodies against protein biomarkers and multiplexed MLs were introduced to recognize the biomarkers through a sandwich protocol.

Results

By analyzing the abundant MLs with laser desorption ionization time-of-flight (LDI-TOF) MS, presence of the biomarkers can be confirmed. As the types of ML increased, a variety of protein biomarkers could be analyzed simultaneously. In this assay, protein biomarkers of epithelial cell adhesion molecule (EpCAM) and cytokeratin 19 (CK19) were selected as exemplified, 2 z mole of which could be specifically and accurately identified. Meanwhile, utilizing this method, EpCAM and CK19 with the amount of 20 z mole were simultaneously detected.

Conclusions

In conclusion, a novel MS based approach has been developed to co-detect multiple biomarkers using numerous MLs-labeled AuNPs. Protein biomarkers, EpCAM and CK19, were simultaneously identified when their amount were lower to 20 z mole. In view of the advantages of ultrahigh sensitivity and specificity, the protocol possesses promising value in cell characterization, cell signaling behavior measurement, and the monitoring of disease formation and progression.

Novel Aspect

By analyzing various MLs with LDI-TOF MS, multiple protein biomarkers could be analyzed simultaneously and sensitively.

References

1. Munge B. S., Stracensky T., Gamez K., DiBiase D., Rusling J. F., *Electroanalysis*, 28, 2644-2658 (2016).
2. Ma D. L., Wang M. D., Liu C. F., Miao X. M., Kang T. S., Leung C. H., *Coord. Chem. Rev.*, 324, 90-105(2016).
3. Ho S. L., Xu D., Wong M. S., Li H. W., *Chem. Sci.*, 7, 2695-2700(2016).
4. Wang L. Y., Rong Q. F., Ma Z. F., *Electroanalysis*, 28, 1692-1699(2016).
5. Du R. J., Zhu L. N., Gan J. R., Wang Y. N., Qiao L., Liu B. H., *Anal. Chem.*, 88, 6767-6772(2016).

TPS-S05 - MS imaging: instrumentation

Chairs: Liam McDonnell, Fulvio Magni

TP-115 / DEVELOPMENT OF OPTICS OF A STIGMATIC IMAGING MASS SPECTROMETER WITH HIGH SPATIAL RESOLUTION

Jun Aoki

Osaka University, Graduate School of Science, Toyonaka

Keywords: imaging mass spectrometry, instrumentation, high spatial resolution

Introduction

Measurement methods of spatial distribution of molecules at cellular-scale are required in many fields. Recently, scanning type imaging mass spectrometry with matrix-assisted laser desorption/ionization is intensively used for biomolecular analysis. However, the spatial resolution of scanning MALDI-IMS is limited to about 10 - 100 μm and inadequate for cellular-scale observation. Therefore, we are developing a stigmatic MALDI imaging mass spectrometer to achieve spatial resolution of sub-micron. [1]

Methods

The experimental apparatus for stigmatic imaging consists of a MALDI ion source with focusing ion optics, TOF mass spectrometer and a time and position sensitive ion detector. Sample molecules ionized by laser irradiation are extracted and accelerated by extraction electrode. Ion distributions at the sample plate are magnified and projected with the ion optical lens system onto the detector.

Results

For a stigmatic MALDI imaging, laser irradiation spot is focused to several hundreds micro-meter diameter with accurate uniformity. This diameter is several times larger than ordinary MALDI, so laser power is required about ten times larger than ordinary one. And for uniform laser irradiation, we develop a laser beam homogenizer system consists of a dynamic rotating diffuser to eliminate laser speckle and a square core optical fiber to form flat-top square profile. We used this system with a high power UV (355 nm) pulse laser (NL202, Ekspla, Lithuania) and evaluated imaging performance by observation of test samples. We observed grid pattern formed by test samples of crystal violet dye covered with fine metal mesh. We verified imaging quality by this experiment and then observed more practical samples of eye tissue of mouse and amoeba cells.

Conclusions

We developed new stigmatic imaging mass spectrometer with a uniform laser beam system. Our evaluation

experiments verified that this apparatus achieved high spatial resolution of 1 micro-meter with a field of view of 500 x 500 micro-meter within observation time of several minutes. This spatial resolution and measurement speed was drastically improved from conventional scanning type imaging mass spectrometers. We also show applications of this apparatus to actual biological samples.

Novel Aspect

We developed optics of stigmatic MALDI imaging mass spectrometer and observed high spatial resolution image.

References

1. Hazama H., Yoshimura H., Aoki J., Nagao H., Toyoda M., et al., J. Biomed. Opt., 16 046007 (2011).

TP-116 / DEVELOPING MS IMAGING IN THE CONTEXT OF CHEMICAL ECOLOGY - LASER ABLATION ELECTROSPRAY IONIZATION ON THREE-DIMENSIONAL SAMPLES

Benjamin Bartels⁽¹⁾ - ***Dirk Hölscher***⁽²⁾ - ***Purva Kulkarni***⁽³⁾ - ***Aleš Svatoš***⁽¹⁾

Research Group Mass Spectrometry/Proteomics, Max Planck Institute for Chemical Ecology, Jena (1) - Ökologische Agrarwissenschaften (OPATS), Universität Kassel, Witzenhausen (2) - Terrestrial Ecology, Netherlands Institute of Ecology (NIOO-KNAW), Wageningen (3)

Keywords: Mass spectrometry imaging, LAESI, profilometry, multimodal, application

Introduction

Laser ablation electrospray ionization (LAESI) is an ionization method applicable for mass spectrometry imaging (MSI) on biological samples [1]. Making use of indigenous water, LAESI is not dependent on externally applied matrices, reducing the sample preparation component of any experiment considerably. By integrating a profilometer, sample topography becomes measurable and employable during laser ablation, allowing access to ecologically relevant samples.

Methods

Various parts of *Gossypium hirsutum* plants were used as a model system to evaluate a custom build LAESI source for MSI [2]. Before conducting MSI experiments, the surface topography was measured and subsequently applied to guide the laser ablation process during the actual MS experiment. Sample preparation was unnecessary, except trimming the sample to fit a microscopy slide. Mass spectra were acquired in positive and negative ion mode on a Synapt HDMS (Waters).

Results

During MSI experiments performed on different ROIs of leaves and stems with our custom-build LAESI source, profilometry was successfully applied to guide the laser ablation. As a result, laser ablation marks were of equal diameter, reproducible and independent of sample surface topography. Specialized metabolites with known insect deterring properties [3], such as gossypol and heliocide isoforms were localized in specialized glands distributed all over *G. hirsutum* plants. By exploiting the measured height profiles, the diameter of laser ablation marks was controllable between 30 µm and ~150 µm in a reproducible manner. Based on experiments on *G. hirsutum*, ecologically relevant interactions between various other plant species and microbial interaction partners were investigated.

Conclusions

Our custom-build LAESI source allows for the inclusion of a profilometry step in MSI experiments. The described set-up is able to deal with challenges posed by ecologically relevant samples of delicate plant-microbial interactions.

Novel Aspect

Analysis of various plant surfaces and plant-microbe interactions without the need for sample preparation or sectioning.

References

1. Nemes P., Vertes A., Analytical Chemistry, 79, 8098-8106 (2007).
2. Bartels B. et al., RSC ADVANCES, 7, 9045-9050 (2017).
3. Bottger G. T., Patana R., J. Econ. Entomol., 59, 1166-1168 (1966).

TP-117 / IMPROVED SENSITIVITY AND SPEED OF ACQUISITION USING QUADRUPOLE BASED MASS SPECTROMETER FOR DESI-MSI

Emmanuelle Claude⁽¹⁾ - Praveen Harapanahalli⁽²⁾ - Emrys Jones⁽³⁾

Waters Corporation, Scientific Operations, Wilmslow (1) - Waters Corporation, Software, Wilmslow (2) - Waters Corporation, Research, Wilmslow (3)

Introduction

DESI-MSI is typically known for the mapping of small molecules directly from tissue sections. DESI MSI has been proven to be successful when used on Time-Of-Flight (TOF) based mass spectrometers in particular for untargeted analysis or mass spectral fingerprinting. Tandem quadrupole mass spectrometers are renowned for their sensitivity and specificity for targeted application using Single Ion Monitoring (SIR) or Multi Reaction Monitoring (MRM) modes of acquisition.

Methods

Experiments were performed on a Xevo TQ-XS (Waters) in SIR or MRM modes of acquisition and a Xevo Q-ToF G2-XS (Waters) in MS or MS/MS modes of acquisition. A 2D DESI stage (Prosolia) DESI was used with spray conditions set 2 μ L/min, 95:5 MeOH:water and a nebulising gas pressure of 5 bar. DESI imaging datasets were subsequently processed and visualised using in-house software tools and High Definition Imaging (HDI) software.

Results

To evaluate the level of sensitivity, dilution series of several drug compounds including Propranolol and Bromofenac were spotted on tissue sections from 0.25ng up to 1 μ g. On the TQ system, MRM mode of acquisition was set with a first transition for the drug and second transition for an endogenous lipid. On the oa-TOF system, MS and MS/MS modes of acquisition were used. The limit of detection for propranolol using the TOF system was 2.5 ng on tissue whereas it was 0.25ng for the DESI TQ. The limit of detection for Bromofenac using the TOF system was 50 ng on tissue in MS and whereas it was 0.25ng on tissue for the DESI TQ. A similar improvement in level detection of drugs was observed from dosed tissue section using the DESI TQ allowing mapping drugs in tissue at therapeutic levels which is not straightforward using TOF mass spectrometer for MSI.

With the improved sensitivity, it was possible to increase the speed of analysis using stage rate from 1 to 30 scans/or the number of MRM, resulting in a decreased dwell time for each pixel.

Conclusions

DESI mounted on TQ mass spectrometer has clearly demonstrated some advantages in an imaging setting such as increased sensitivity allowing drug compounds to be mapped at therapeutic level from dosed tissue sections.

Novel Aspect

Mapping of drug compounds at therapeutic level using highly sensitive and selective DESI mounted on triple quadrupole mass spectrometer.

TP-118 / ATMOSPHERIC PRESSURE MALDI IMAGING ORBITRAP MS USING A MASSTECH APMALDI UHR ION SOURCE

Gilles Frache - *Dana El Assad*

Luxembourg Institute of Science and Technology (LIST), Materials Research and Technology, Belvaux

Keywords: targeted and untargeted AP-MALDI imaging

Introduction

The interest of scientist for MALDI imaging is growing since the last decade and latest developments has led to a significant impact in the pharmaceutical and cosmetics field, but also in materials research. MALDI imaging is now enabled on high resolution Orbitrap instruments by means of Atmospheric Pressure MALDI sources available on the market.

Methods

An AP-MALDI UHR ion source (Masstech Inc.) is coupled to a LTQ/Orbitrap Elite high resolution mass spectrometer (Thermo Scientific) was used for the development of targeted and untargeted imaging Mass Spectrometry experiments. Samples were coated with matrices using a HTX-TM sprayer.

Results

Several modes of operation of the AP-MALDI Orbitrap setup are presented based on samples from various types and sizes for the detection of targeted and untargeted molecules. The instrument allows for both large and small scale imaging (using Pixel map or Constant speed raster modes).

The information extracted from each image pixel may be adjusted to either untargeted (full scan HRMS) or targeted (MS/MS or MSn). These various experiments may performed simultaneously in both analyzers in each location of the sample.

Conclusions

We show the characteristics and applications of a Masstech AP-MALDI UHR ion source to enable targeted and untargeted MALDI imaging capabilities down to 10 micron lateral resolution, when associated to a LTQ/Orbitrap mass spectrometer. This system allows to rapidly switch from the MALDI configuration to LCMS configuration.

Novel Aspect

Easily interchangeable AP-MALDI (Masstech UHR ion source) for targeted and untargeted MALDI HR-MS or –MSn imaging capabilities down to 10 micron lateral resolution.

TP-119 / HYBRID SIMS: A NEW SIMS INSTRUMENT FOR HIGH RESOLUTION ORGANIC IMAGING WITH HIGHEST MASS-RESOLVING POWER AND MS/MS

Sven Kayser⁽¹⁾ - *Alexander Pirkl*⁽²⁾ - *Rudolf Moellers*⁽²⁾ - *David Scurr*⁽³⁾ - *Nichola Starr*⁽³⁾ - *Ewald Niehuis*⁽⁴⁾

IONTOF GmbH, Sales, Muenster (1) - IONTOF Technologies GmbH, R&D, Muenster (2) - University of Nottingham, School of Pharmacy, Nottingham (3) - IONTOF Technologies GmbH, Management, Muenster (4)

Keywords: SIMS, Imaging, MS/MS, Orbitrap

Introduction

SIMS offers the possibility to acquire chemical information from submicron regions. This capability has been especially intriguing for researchers with life science applications. The vision to image molecules has

been driving instrumental development. While new ion sources expanded the usability, SIMS analysers lacked the required mass resolution, mass accuracy and MS/MS capabilities.

Methods

To specifically address the requirements in the life science field a powerful new Hybrid SIMS instrument [1] was developed in a research project by IONTOF and Thermo Fisher Scientific™, in close cooperation with the National Centre of Excellence in Mass Spectrometry Imaging, GlaxoSmithKline, and the University of Nottingham, following Prof Ian S Gilmore's original idea.

Results

In this contribution we will demonstrate the possibilities of the new instrumental design, which combines an Orbitrap™-based Thermo Scientific™ Q Exactive™ HF mass analyser with a high-end ToF-SIMS system from IONTOF GmbH. First application data of high resolution SIMS spectrometry, MS/MS analyses, high resolution imaging of tissues and depth profiles of biological samples will be presented. Especially the power of combining the two mass analysers on a single instrument – time-of-flight analyser for fast images with highest spatial resolution and orbital trapping analyser for highest mass resolution and mass accuracy – will be highlighted. results also cover single beam depth profiles through porcine skin samples, which allow for a clear skin layer elucidation following different molecular ion signals.

Conclusions

The new Hybrid SIMS instrument combines for the first time highest mass resolution (> 240,000) and highest mass accuracy (< 1 ppm) and full MS/MS capabilities with high resolution cluster SIMS imaging, making SIMS an even more viable technique for the analysis of life science samples.

Novel Aspect

The combination of the imaging capabilities of a TOF analyser with the unique performance of the Orbitrap™ provides a new level of SIMS information.

References

1. Passarelli et al., Nature Methods, 14, 1175–1183 (2017)

TP-120 / DESORPTION ELECTROSPRAY IONIZATION TRIPLE QUADRUPOLE MASS SPECTROMETRY: TARGETED DRUG AND METABOLITE IMAGING IN BIOANALYTICAL RESEARCH

Lieke Lamont⁽¹⁾ - **Gert Eijkel**⁽¹⁾ - **Emrys Jones**⁽²⁾ - **Bryn Flinders**⁽¹⁾ - **Shane Ellis**⁽¹⁾ - **Tiffany Porta**⁽¹⁾ - **Rob J. Vreeken**⁽³⁾ - **Ron M. A. Heeren**⁽¹⁾

Maastricht University, Maastricht Multimodal Molecular Imaging (M4I) Institute, Maastricht (1) - Waters Corporation, Waters Corporation, Wilmslow (2) - Janssen Pharmaceutica, Janssen R&D, Beerse (3)

Keywords: Targeted MS imaging, multiple reaction monitoring, neutral loss scanning, precursor ion scanning, drug discovery

Introduction

Drug discovery aims at finding suitable drug candidates that are being metabolized and distributed to the target site, bind to the biological receptors, and block or stimulate a metabolic reaction. [1]Mass spectrometry imaging has established itself as a powerful tool to provide spatial information on drugs and metabolites but suffers from low sensitivity and irreproducibility due to background ions from biological sample matrix. [2]

Methods

Our work reports a novel targeted MS imaging approach using desorption electrospray ionization (DESI)

combined with multiple reaction monitoring (MRM), neutral loss (NL), precursor ion (PREC), and product ion (PROD) scans for the analysis of drugs and lipids in pre-clinical studies. MRM uses a precursor ion-product ion transition to provide the selectivity required to separate isobaric species based on different product ions.

Results

For this purpose, a set of experiments is executed on a pre-clinical drug study in which multiple dog liver tissues were analyzed using DESI prior to PROD, NL, and MRM modes to target, map and identify two potential drug candidates. The accumulation of the drugs correlated with histological lesions (already annotated by the pathologist) indicates severe fibrosis which resulted in stopping these compounds from further development. In addition, we further assessed the potential of MRM, NL, and PREC imaging modes by monitoring several endogenous lipid species in consecutive tissue sections, to complement the investigation of the toxic effect. The improved selectivity and sensitivity of this technology was investigated by the comparison of NL and PREC scanning modes with the MRM mode of low abundant lipid species in rat brain.

Conclusions

Our approach opens doors for the pharmaceutical industry to target and image drugs and metabolites with high selectivity without the need for expensive high mass resolution imaging technologies. The use of a triple quadrupole MS not only offers higher selectivity and sensitivity but also allows to explore a (semi) quantitative strategy due to its high dynamic range.

Novel Aspect

The selective capability of a triple quadrupole is now combined with mass spectrometry imaging for the first time by using targeted NL, PREC, and MRM imaging.

References

1. Prideaux B., Stoeckli M., Journal of Proteomics, 75, 4999-5013 (2012).
2. Laskin J., Lanekoff I., Analytical Chemistry, 88, 52-73 (2016).

TP-121 / OPEN DEVELOPMENT KIT FOR BUILDING AN AMBIENT MASS IMAGING ROBOT WITH MICROMETRIC RESOLUTION

Abigail Moreno-Pedraza⁽¹⁾ - **Cesaré Ovando-Vázquez**⁽²⁾ - **Héctor Guillén-Alonso**⁽¹⁾ - **Ignacio Rosas-Román**⁽¹⁾ - **Sebastian Gibb**⁽³⁾ - **Robert Winkler**⁽¹⁾

CINVESTAV, Biotechnology and Biochemistry, Irapuato (1) - IPICYT, National Supercomputing Center, San Luis Potosí (2) - Department of Anaesthesiology and Intensive Care, University Medicine Greifswald, Greifswald (3)

Keywords: ambient ionization, ambient imaging, 3D printing, low temperature plasma, instrumentation

Introduction

Ambient ionization has revolutionized mass spectrometry, because samples can be analyzed under atmospheric condition without prior preparation, including imaging [1]. Technologies such as 3D printing help to easily share, implement and reproduce ambient ionization probes [2]. Inspired by this technology, we decided to generate an automated sampling platform based on RepRap systems for high throughput analyses and imaging to further enhance the utility of AIMS.

Methods

RepRap technology was used for the ambient sampling platform. The system was designed in our

laboratory and built by MakerMex (Gto. Mexico). Additional accessories to adjust ionization sources and optical tools were designed using Blender and printed using a MM1 Beta RepRap printer with an acrylonitrilebutadiene styrene (ABS) polymer. For system control and data analysis, the software rmsiGUI was built based on the free statistics software R using:shiny package.

Results

We developed a complete kit for an ambient sampling platform, consisting of -robot -software -adapters. The sampling platform can incorporate new tools/accessories depending on the applications. The system can be used for imaging and high throughput analyses. The theoretical resolution of our platform is 12.5 μm . However, the natural diffusion of the compounds generates a practical resolution of 200 μm . We used a PlasmaChip (NovionX, Germany) to automate sample analyses and for ambient imaging we used a 3D-plasma probe alone or in combination with laser desorption (LD). The rmsiGUI software has tabs for different functions, generates a G-code to control the robot movements, creates imzML files for mass spectrometry imaging, uploads a photo from a sample and marks a region of interest (ROI). Data can be analyzed by multiple statistical tools to identify metabolic markers.

Conclusions

We demonstrated the functionality of an open, modular ambient sampling platform for automated sampling and mass imaging. The rmsiGUI is able to perform a complete analysis from MSI data. The design of the sampling robot and the rmsiGUI software are freely available for non-commercial uses.

Novel Aspect

We present an open platform for analysis under ambient conditions based on RepRap technology which enables the user to design new features based on his specific needs.

References

1. Domin M., Cody R., Royal Society of Chemistry, 508 (2015).
2. Martínez-Jarquín, S., Moreno-Pedraza, A., Guillén-Alonso, H., Winkler, R., Anal. Chem., 88, 6976–6980 (2016)

TP-122 / BIG DATA POST-PROCESSING FOR HIGHER THROUGHPUT MALDI/DESI FTMS IMAGING

Konstantin O. Nagornov⁽¹⁾ - ***Anton N. Kozhinov***⁽¹⁾ - ***Pieter C. Kooijman***⁽²⁾ - ***Edith Nicol***⁽³⁾ - ***David P. A. Kilgour***⁽⁴⁾ - ***Shane Ellis***⁽²⁾ - ***Ron M. A. Heeren***⁽²⁾ - ***Yury O. Tsybin***⁽¹⁾

Spectroswiss, Biomolecular Mass Spectrometry, Lausanne (1) - Maastricht University, Imaging Mass Spectrometry, Maastricht (2) - Ecole Polytechnique, Laboratoire de Chimie Moléculaire, Palaiseau (3) - Nottingham Trent University, Chemistry and Forensics, Nottingham (4)

Keywords: FTMS, MALDI/DESI, imaging, transient, software

Introduction

High mass resolution FTMS MALDI and DESI imaging techniques are widely employed for analyzing biological tissues, plants, and oils. Real-life sample complexity, tiny quantities of analytes, and long experimental times require increased sensitivity, resolution, and throughput. Here, we describe innovations in FTMS transient acquisition and processing algorithms and their integration into big data analysis workflows for leveraging FTMS imaging applications.

Methods

MALDI and DESI imaging data were acquired with commercial ICR and Orbitrap FTMS instruments. Mass spectra and transients were recorded in parallel using original instrument manufacturer and in-house developed data acquisition hardware and software. Data post-processing was performed using custom

software tools and algorithms, including absorption mode FT (aFT) and least-squares fitting (LSF), using multicore and graphic card processing capabilities.

Results

We evaluated the key benefits resulting from the use of a novel big data processing workflow. Disabling on-the-fly FTMS data processing (including FT, peak picking, and recording data on a drive) reduced the total time of measurements up to 2-fold per image. Additionally, the scan rate was further increased via an optimized use of overhead times and application of aFT (allows approx. half the transient time compared to magnitude mode processing) or LSF for targeted imaging. The LSF demonstrated superior performance to magnitude and aFT in terms of required minimum transient length. LSF processing, being implemented on graphic cards, reduced data processing times to practically-useful acceptable levels. Increases in sensitivity and identification confidence were achieved by transient or full-profile mass spectra summation over a region of interest and summation across multiple experiments, *e.g.* 3D MALDI/DESI. Mass accuracy was improved using an internal pixel-by-pixel non-linear mass re-calibration based on binomial averaging.

Conclusions

Increasing the rate at which FTMS transients are generated, while storing them for post-processing and re-analysis, appears attractive for higher-throughput MALDI/DESI imaging applications. The aFT and LSF methods allow the use of shorter transients, without compromised accuracy in identification of molecules (*e.g.*, lipids). Increased sensitivity enhances the confidence and the number of lipids identified and quantified per MALDI/DESI imaging experiment.

Novel Aspect

Advanced algorithms of time-domain signal (transient) acquisition and post-processing, and big-data analysis boost FTMS performance for MALDI/DESI imaging applications.

TP-123 / PATERNÒ-BÜCHI PHOTO-DERIVATIZATION PROTOCOL FOR MALDI-MS IMAGING OF CARBON-CARBON DOUBLE BOND POSITIONAL ISOMERS OF PHOSPHO- AND GLYCOLIPIDS IN BIOLOGICAL TISSUES

Antonín Bednařík - Stefan Bölsker - Jens Soltwisch - Klaus Dreisewerd

University of Münster, Institute for Hygiene - Biomedical Mass Spectrometry, Münster, Germany

Keywords: Paternò-Büchi reaction, Benzaldehyde, MALDI-MS imaging, Double bond positional isomers, MALDI-2

Introduction

Current MALDI-MS imaging (MALDI-MSI) techniques visualize the global distribution of phospho- and glycolipids in tissue sections, but are generally blind to the physiologically relevant carbon-carbon double bond (C=C) position. Here we present a fully MALDI-(2)-MSI compatible “on-tissue” Paternò-Büchi (PB) photo-derivatization protocol that enables mapping the distribution of C=C positional isomers of numerous phospho- and glycolipids in mouse brain.

Methods

Benzaldehyde (BA; *m* monoisotopic 106.04 Da) vapor was introduced into a reaction chamber and condensed on the sample (sections of murine brain) at ~6°C. PB reaction was triggered for 1-2 min by use of a 254 nm-UV lamp. Subsequently, BA was allowed to evaporate at ambient temperature. After coating with matrix, samples were analyzed at 25 μm pixel size with a Q Exactive plus Orbitrap that was equipped with a dual MALDI/ESI injector and a postionization laser for MALDI-2 [1,2].

Results

Making use of our combined PB-MALDI-2 approach, numerous phospho- and glycolipids produced clear

$[M+BA+H]^+$ ion signals in the positive ion mode. Overall PB photo-derivatization yields were in the 5–30% range (depending on lipid class). Low-energy CID of the protonated PB products generated clear MS/MS spectra and characteristic ion pairs exhibiting a mass difference of 74.05 Da. These unambiguously revealed the C=C positions for numerous singly and poly-unsaturated phospho- and glycosphingolipids in mouse cerebellum. Several C=C positional isomers of phosphatidylcholine and, for the first time, phosphatidylserine species were found to be expressed with highly differential levels in the white and gray matter areas of cerebellum. In contrast, C=C isomers of galactosylceramides were found with equal distributions in the myelin-rich white matter. Importantly, at the applied lateral resolution of 25 μm , no tissue damage or analyte delocalization, which could be associated with our PB protocol, was observed.

Conclusions

Our PB-MALDI-2-MS/MS protocol provides a valuable glimpse into phospholipid organization complexity in tissues with currently unknown biological background. Compared to the also just introduced ozon-based derivatization protocol (OzID) for differentiation of C=C positional isomers [3], our approach has the advantage to widen the class of amenable lipids; moreover, use of any standard MALDI-MSI instrument (also without postionization module) is possible.

Novel Aspect

Paternò-Büchi reaction-based protocol to visualize the distribution of C=C positional isomers of phospho- and glycolipids by MALDI-MS imaging at high lateral resolution

References

1. Soltwisch J., Ketting H., Vens-Cappell S., Wiegelmann M., Müthing J., Dreisewerd K., *Science*, 348, 211–215 (2015)
2. Ellis S. R., Soltwisch J., Paine M. R. L., Dreisewerd K., Heeren R. M. A., *Chemical Communications*, 53, 7246–7249 (2017)
3. Paine M. R. L., Poad B. L. J., Eijkel G. B., Marshall D. L., Blanksby S. J., Heeren R. M. A., Ellis S. R. *Angewandte Chemie International Edition*, (2018) <https://doi.org/10.1002/ange.201802937>

TP-124 / HEAT-ASSISTED LASER ABLATION ELECTROSPRAY IONIZATION (HA-LAESI) FOR MS IMAGING OF ANIMAL TISSUE

Juha-Pekka Hieta - Riikka-Marjaana Räsänen - Jaakko Kopra - Heikki Rääkkönen - Risto Kostainen - Tiina Kauppila

University of Helsinki, Faculty of Pharmacy, Helsinki, Finland

Keywords: HA-LAESI, LAESI, imaging, tissue, laser

Introduction

Heat-assisted laser ablation electrospray ionization (HA-LAESI) [1] and laser ablation electrospray ionization (LAESI) [2] are ambient MS techniques for the sensitive analysis of compounds with low molecular weight (<1000) from water-rich samples. Previously, HA-LAESI has been showed to be more efficient in ionizing nonpolar and neutral compounds than LAESI, and thus HA-LAESI is expected to show complementary data in the imaging of animal tissue.

Methods

Imaging measurements were carried out using a Bruker micrOTOF MS. In HA-LAESI measurements, tissue was ablated using an OPOTEK IR Opolette 2940 laser, and a heated capillary extension was used to add heat to electrospray which was emitted with a typical 100 μm I. D. stainless steel needle. The tissue sections were cooled to -4°C using a Peltier element and they were rastered inside a N_2 filled box. Waters Xevo QTOF with a ACQUITY UPLC was used for compound identification.

Results

Previously, a simple technique for improving the focus of OPO IR lasers has been demonstrated by solely increasing the laser beam travel distance [3]. The new focusing technique was combined with HA-LAESI to show heat maps of tissue sections with <50 μm pixel size. The HA-LAESI imaging technique was applied to several types of whole tissues, such as mouse brain, heart and kidney. Some of the detected analytes were identified using Waters' UPLC/MS system. E. g. from mouse brain tissue, HA-LAESI showed spectra covered with high intensity phospholipid, cholesterol, ceramide, and cerebroside signals.

Conclusions

HA-LAESI was found to be a suitable technique for high resolution imaging of animal tissues, and it showed to be more efficient in ionizing nonpolar compounds than LAESI.

Novel Aspect

High resolution HA-LAESI imaging of animal tissues

References

1. Vaikkinen A. et al., Anal. Chem., 85, 177–184 (2012).
2. Nemes P. et al., Anal. Chem., 79, 8098–8106 (2007).
3. Hieta JP. et al., J. Am. Soc. Mass Spectrom., 28, 1060–1065(2017).

TP-125 / PERFORMANCE-ENHANCED DATA ACQUISITION UTILIZING REAL-TIME DSP

Ulrik Lindblad - Carsten Watolla

Teledyne Signal Processing Devices Sweden AB, N/A, Linkoping, Sweden

Keywords: data acquisition, digitizer, digital signal processing, pattern noise reduction, baseline drift

Introduction

Accurate pulse detection is vital in data acquisition systems for mass spectrometry but problems such as component aging and temperature variations typically reduce the accuracy of the detection and may generate misleading results. Drift in baseline level and pattern noise from interleaved analog-to-digital converters (ADCs) may result in missed pulses or incorrect detection of false pulses.

Methods

Data acquisition modules from Teledyne SP Devices utilize on-board real-time signal processing to remove baseline drift and reduce ADC pattern noise. The Digital Baseline Stabilizer (DBS) operate in the background and perform highly accurate baseline adjustment to a user-defined target value. ADC pattern noise reduction is achieved by continuous background calibration that corrects for gain, offset and time skew mismatch errors.

Results

Baseline adjustment and pattern noise reduction results in significantly improved detection accuracy. Continuous background correction allows the user-controlled trigger level to be valid over time without any need for manual calibration procedures. Baseline drift and pattern noise can be significant sources of noise – the latter with mismatch errors ranging up to a hundred ADC codes [1]. DBS reduces the pattern noise to a level corresponding to 22 bits resolution.

Conclusions

High-precision pulse detection requires estimation and correction of imperfections and mismatch errors in analog electronic components. Baseline drift in data acquisition for mass spectrometry can lead to misleading pulse characteristics and peak values as they are determined against an incorrect baseline level and ADC pattern noise can further significantly contribute to false readings.

Novel Aspect

DBS shows how expertise in signal processing, hardware knowledge, and application know-how is combined to create an optimized system solution.

References

1. Teledyne SP Devices, The art of pulse detection (2018).

TP-126 / ELECTROSPRAY IONIZATION/ATMOSPHERIC PRESSURE INTERFACES: COMPREHENSIVE SIMULATION OF GAS FLOW, SPACE CHARGE AND ELECTRIC FIELDS

Philipp Krahl⁽¹⁾ - *Laurent Bernier*⁽¹⁾ - *Stephan Rauschenbach*⁽²⁾ - *Julius Reiss*⁽¹⁾

Technical University Berlin, ISTA, Berlin, Germany (1) - University of Oxford, Dept. of Chemistry, Oxford, United Kingdom (2)

Keywords: ESI, simulation, gas flow, space charge, ion funnel

Introduction

In electrospray ionization (ESI) ions are often transferred through a narrow capillary from ambient to a lower pressure stage, where complex ion optics like ion funnels are used to form the ion beam. The particle transport is governed by electric forces, space charge and by the gas flow. Their interplay determines the behavior and in particular the efficiency of the source. However, most simulations include only some of these effects, often neglecting the gas dynamics.

Methods

We simulate the gas flow of the capillary and the first pumping stage including an ion funnel. The resulting flow fields are used to simulate the ion transport by a self-written software. It allows to treat the space charge effects and collimation by the rf- and dc-fields of the ion funnel by first principles and to quantify the different effects. It also reveals the conditions the ions experience, like temperature, pressure and electric fields, during transfer.

Results

We find a strong influence of the gas flow on the ion transmission. High transmission rates recently observed can be reproduced by the simulation and understood, with numerical results comparing quantitatively to the experiment. The key factor to obtain high transmission in a transfer capillary is transport in a laminar flow. High transmission is further obtained for low net currents (low space charge repulsion) and good pre-collimation of an electrospray before entering the transfer capillary, as for instance achieved by recently presented funnel inlets. But also under ideal laminar conditions, space charge effects determine the maximal current transmitted through the capillary. Time resolved simulations with many individual particles allow an in-depth description especially in the ion funnel. When designing an ESI sources all the different aspects need to be taken into account.

Conclusions

Gas flows in the high pressure part of an ESI interface are important for ion transport and deserve higher attention due to their large potential for optimizing ion sources. Space charge effects become important for high currents and can partly be mitigated by an improved extraction and collimation at the inlet by the gas flows. Better understanding of the interplay of space charge and gas flow is the key to high current and high transmission ion sources.

Novel Aspect

A comprehensive simulation of the transfer capillary and the first stage with an ion funnel reveals detailed insight into important aspects of the ion transport.

TPS-S06 - Post Translational Modifications

Chairs: Joseph A. Loo, Fulvio Magni

TP-127 / IMPROVEMENT OF PROTEOMIC ANALYSIS TARGETING MULTIPLE POST-TRANSLATIONAL MODIFICATIONS USING IMMUNOPRECIPITATION WITH PAN-PTM ANTIBODIES

Yuichi Abe - Asa Tada - Junko Isoyama - Takeshi Tomonaga - Jun Adachi

National Institute of Biomedical Innovation, Health, and Nutrition, Proteome Research Project, Osaka-fu

Keywords: post-translation modification, systems oncology, immunoprecipitation, PTM crosstalk

Introduction

It is widely known that post-translational modification (PTM) is important for regulation of protein activity and several cellular functions. Therefore, abnormality in those modification is closely related to oncogenesis. Thus, understanding of global PTM status significantly contribute to cancer biology and drug discovery against cancer. So far, several PTM proteomics were conducted by using immunoprecipitation with pan-PTM antibodies.

Methods

One of examples is phosphotyrosine proteomics. When performing large-scale phosphoproteomics, phosphorylated tyrosine sites account for less than 1%. To increase identification of phosphotyrosine sites, we optimized enrichment step of phosphotyrosine peptides with immunoprecipitation, and achieved highly sensitive detection of phosphotyrosine peptides (more than 1,000 sites) from samples without pervanadate treatment [1].

Results

In this study, we examined optimization for enrichment step of peptides containing other PTM sites. Newly developed protocol is commonly applicable for enrichment of several types of PTMs. We successfully measured several PTM status from single sample.

Conclusions

We concluded that optimization for enrichment of each type of PTM peptides is crucial for comprehensive proteomics of multiple PTM.

Novel Aspect

This technique would lead to understand complexed crosstalk between multiple PTM signalings and find novel therapeutic targets.

References

1. Abe Y et al., Journal of Proteome Research, 16, 1077-1086 (2017).

TP-128 / THE ROLE OF THIOLS IN REDOX BIOLOGY: HPLC-MS CHARACTERIZATION OF CYSTEINE DERIVATIVES GENERATED BY INTERACTION WITH COLD PLASMA JETS, POWERFUL SOURCES OF NITROGEN AND OXYGEN SPECIES

Giuliana Bruno - Jan-Wilm Lackmann - Sebastian Wenske - Sander Bekeschus - Klaus-Dieter Weltmann - Thomas von Woedtke - Kristian Wende

Leibniz Institute for Plasma Science and Technology (INP Greifswald), ZIK plasmatis, Greifswald

Keywords: thiols, redox biology, cold physical plasmas

Introduction

The proteome is a sensor and mirror of the redox state in mammalian cells [1]. Cysteine is the amino acid most sensitive to oxidation since its thiol group forms different biologically relevant chemotypes [2]. In this study, possible derivatives of cysteine formed by impact with cold physical plasma as means to produce oxygen and nitrogen centered reactive species (RONS) [3] were surveyed by high resolution MS/MS in order to predict the potential impact of plasma in vivo.

Methods

Different concentrations of cysteine standard solutions were treated by various cold plasma jets with respective parameter variation. The unknown compounds were then separated using an HILIC column in an HPLC system, and identified by using a high resolution q-TOF mass spectrometer. A relative quantification was made by relating the intensities in each spectra to the one of the L-cysteine-¹³C₃, ¹⁵N at known concentration.

Results

A large number of oxidation products from cysteine were identified and the potential pathway of formation relative to RONS generation by the plasma was described. Modifications were mostly oxygen driven and located on the thiol moiety, with the cysteine sulfonic acid as major, but not ultimate, compound. By distinct modulation of the plasma parameters, *e.g.* nitrogen shielding of the effluent or addition of water to the feed gas, the formation of specific cysteine bioactive derivatives such as S-nitrosocysteine, cysteine-S-sulfonate, or cystine di- and tetra-sulfoxides were increased. By varying the plasma parameters, cysteine's biochemical impact could be modulated ultimately reflecting back on the biological impact of cold plasma treatment in redox related conditions, *e.g.* wound healing[3].

Conclusions

Cold plasmas are valid tools to produce RONS and simulate molecular changes occurring in oxidative environments. It was observed that cysteine product profile reflects treatment conditions, allowing cysteine variants with potential relevance in redox signaling to be characterized. These observations imply a profound impact on protein thiol groups when treated by oxidizing conditions like plasma in vivo, with consequences on cellular redox signaling pathways.

Novel Aspect

Rare but potentially biologically relevant cysteine derivatives were observed after physical plasma treatment, allowing a forecast of thiol oxidation products in vitro.

References

1. Y. M. Go, D. P. Jones, Redox Biol. 2, 358-360 (2014).
2. N. O. Baez, J. A. Reisz, C. M. Furdul. Free Radical Biol. Med., 80, 191-211 (2016).
3. K. D. Weltmann, T. von Woedtke. Plasma Phys. Controlled Fusion, 59, 014031 (2017).

TP-129 / A RAPID, MULTIPLEXED KINASE ACTIVITY ASSAY USING ITRAQ LABELING AND MALDI-TOF/TOF MS

Chao-Jung Chen - Yu-Ching Liu

China Medical University, Taichung

Keywords: multiplexed, kinase assay, MALDI-TOF/TOF, insulin stimulation, inhibitor

Introduction

Substrates from chemically synthesized peptides have been used for in vitro phosphorylation using purified

kinase or cell lysates. Because of its high specificity and high throughput, MALDI-TOF-MS has been a powerful platform for the detection of phosphorylated peptide substrate from enzymatic reactions. However, for screen assay, a direct readout and comparison among different experimental conditions without using internal standard would be preferred.

Methods

The peptide substrates labeled with 8-plex iTRAQ reagents were separately added in 8 different experimental conditions. For kinase activity comparison among 8 samples, the lysates of each experimental conditions containing individual iTRAQ-labeled peptide substrates were pooled for protein removal using SPE, followed by enrichment of phosphorylated peptide substrates using TiO₂ plate and MALDI-TOF analysis.

Results

We found that abundant proteins in cell lysate can severely affect the purification efficiency of phosphorylated peptides on TiO₂ plate. Therefore, SPE approach was first optimized to remove abundant proteins without significant influence of the recovery of phosphorylated peptides from the TiO₂ plate. To enable our platform applicable in insulin stimulation and cancer drug inhibition, a list of the peptide substrates was designed. The protein amount and incubation time for kinase assay of insulin stimulation and drug inhibition were optimized. By labeling peptide substrates with 8-plex iTRAQ reagents, protein kinase activities in 8 samples can be directly compared by the mass tags on their fragmented ion spectra. The kinase activity results from our multiplexed-kinase activity MALDI-TOF assay are consistent with the results from western blots analysis.

Conclusions

We have successfully developed a rapid and quantitative measurement of multikinase activity based on MALDI-TOF. In this assay we combined 8-plex iTRAQ labeled peptide substrates for quantification with a phosphorylated peptide substrate purification plate to rapidly determine multikinases activities in cell lysates. The effect of insulin stimulation and an inhibitory drug on the cellular protein kinase activity were successfully evaluated.

Novel Aspect

This study is the first report to use iTRAQ-labeling strategy combined with TiO₂ plate and MALDI-TOF/TOF for sensitive and high throughput multikinase activity analysis.

TP-130 / USING SUB-RANKED DATABASE MATCHING SCORES FOR IMPROVING THE PEPTIDE AND PROTEIN IDENTIFICATION PERFORMANCE

Yet-Ran Chen - Ying-Lan Chen - Wei-Hung Chang

Academia Sinica, Agricultural Biotechnology Research Center, Taipei

Keywords: Peptidomics, Database Search, Proteomics

Introduction

In peptidomics and proteomic study, a large portion of the MS/MS spectra acquired by the data dependent acquisition (DDA) approach cannot confidently match to a specific peptide sequence. To improve the peptide and protein identification performance, a new database search approach was proposed. Performance and application of this approach were studied and discussed.

Methods

The processed DDA analysis spectra were searched against a protein database using semi-tryptic peptide specific (sTPS) and non-enzyme specific (NES) setting in the identification tryptic peptides and endogenous

peptide, respectively. The 1st to the 10th peptide matching scores for each of the MS/MS spectrum were used to develop new peptide matching score.

Results

In this study, we showed that the score distributions of 1st and 2nd ranked hits were similar if the 1st ranked hits matched to random sequences but distributions were differentiated if the 1st ranked hits matched to forward sequences. Because the differences of the 1st and 2nd ranked score for random matched hits are smaller than the true positive hits, the use of delta score (DS) can make the true positive hits to be more discriminated from the random hits. We further developed a new scoring method named confidence scoring (CS), which considered the 1-10th matching scores of each spectrum. For the identification of tryptic peptide using NES search, the peptide identification sensitivity of CS was 12% and 2% higher than the Mascot score and DS, respectively. In the sTPS search, the number of the identified peptides of CS and DS was 12% and 25% higher than the Mascot score, respectively.

Conclusions

A new approach for performing the tryptic and endogenous peptide identification was developed. This new approach not only can improve the identification performance of endogenous peptides but also on tryptic peptides. In the analysis Tomato peptidomics, we demonstrated that the majority of the endogenous plant peptides are derived from the C-terminal of precursor proteins.

Novel Aspect

This method does not require complex computational steps and can significantly improve the sensitivity for detecting tryptic and endogenous peptide.

TP-131 / SEPARATION OF COMBINATORIAL PTMS ON HISTONE H3 TAIL WITH AN ADVANCED DATA ANALYSIS AND ION MOBILITY

Yoko Chikaoka⁽¹⁾ - *Hiroki Ueda*⁽²⁾ - *Ryo Kajita*⁽³⁾ - *Kazuki Yamamoto*⁽¹⁾ - *Lee Sujin*⁽¹⁾ - *Toshiya Tanaka*⁽⁴⁾ - *Takeshi Kawamura*⁽¹⁾

Isotope Science Center, The University of Tokyo, Tokyo, Japan (1) - Fujitsu Limited, The University of Tokyo, Biological Data Science Division, RCAST, Tokyo (2) - Bruker Japan K. K., Yokohama (3) - The University of Tokyo, Laboratories for Systems Biology, RCAST, Tokyo (4)

Keywords: protein modifications, ion mobility, data analysis, histone, middle-down

Introduction

Mass spectrometry is suitable to investigate protein post-translational modifications (PTMs). Any modifications on any proteins are detectable as mass shifts. Combinations of the PTMs are often called as 'codes' and they are thought to have different biological meanings [1]. Histones are representative proteins that have the code. Since the codes can be formed between PTMs on distant amino acids, typical LC-MS/MS analysis using trypsin often loses their combinatorial information.

Methods

Histone H3 was purified from cells by histone extraction and HPLC fractionation. Histone H3 tail 1-49 amino acid long peptides are digested with AspN. LC-MS/MS with electron transfer dissociation (ETD) was performed using the Orbitrap Fusion with high-resolution MS/MS and Ion mobility separation-MS/MS analysis was performed using the timsTOF Pro. Data was analyzed with MASCOT and a developing software in house.

Results

Histone H3 tails have multiple positional isomers of PTMs that is difficult to separate by liquid chromatography. Consequently, the obtained MS/MS spectra had the mixed ones of the positional isomers

of modifications. An in-house software has been developed by using de Bruijn Graph to resolve the mixed spectra.

We are examining the effects of ion mobility separation [2]. We have obtained several species that have different drift time from one precursor, which might be positional isomers.

Conclusions

Advanced data analysis and ion mobility separation might improve combinatorial PTM identification.

Novel Aspect

We are developing novel software to improve PTM identification of histone middle-down analysis.

References

1. Aebersold R. et al., Nat Chem Biol., 14, 206-214 (2018)
2. Meier F. et al., J. Proteome Res., 14, 5378-5387 (2015)

TP-132 / QUANTITATIVE PHOSPHOPROTEOMICS IDENTIFIES ALTERATIONS OF DOUBLECORTIN-LIKE KINASE 1 (DCLK1) IN NEURAL PROGENITOR CELL DIFFERENTIATION

Alireza Dehghani - Markus Goedderz - Volkmar Gieselmann - Dominic Winter

Institute of Biochemistry and Molecular Biology, University of Bonn, Bonn

Keywords: Quantitative phosphoproteomics, Neurosphere differentiation, LC-MS/MS, DCLK1, BioID

Introduction

Oligodendrocytes are the myelinating cells of the central nervous system and originate from neural stem cells through differentiation mechanism which is only incompletely understood. [1] It has been shown that kinases play important roles in this differentiation process. [2-3] Therefore, in this study, we performed a large-scale phosphoproteomics experiment to investigate phosphoproteome alterations during the differentiation process.

Methods

Samples were harvested at 3 time points of differentiation, labeled using 3-plex dimethyl, and separated into 12 fractions using SCX chromatography. Phosphopeptides were enriched and analyzed by LC-MS/MS. For follow-up experiments, proximity biotinylation by DCLK1-BioID fusion proteins was performed: 8 BioID constructs were generated resembling distinct splice variants of DCLK1. Samples were biotin- and phosphopeptide-enriched and analyzed by LC-MS/MS.

Results

Neural progenitor cells were differentiated over 2 weeks to oligodendrocyte precursor cells. Samples were harvested at three time points: 0, 7, and 14 days of differentiation. Harvested samples were lysed, digested, dimethyl labeled and fractionated using LC-SCX fractionation. Phosphopeptides were enriched using TiO₂ and analyzed by LC-MS/MS. 6963 unique phosphopeptides were quantified covering 9053 phosphorylation sites. In total, 1118 phosphosites were up- or downregulated at least 2-fold. Among the most highly upregulated candidates, several phosphorylated peptides of Doublecortin Like Kinase 1 (DCLK1) were found. Phosphomimetic studies showed that DCLK1 phosphorylation alters its stability. We therefore employed BioID [4]-DCLK1 fusion proteins to investigate DCLK1 substrates and interaction partners. In these experiments the phospho-enriched fractions revealed 2574 phosphosites representing 2086 phosphopeptides from 1229 phosphoproteins while the flow through yielded 33623 peptides from 4616 proteins.

Conclusions

9 residues of DCLK1 were found to be differentially phosphorylated during the differentiation process. Site

directed mutagenesis revealed that phosphorylation in the SP-rich domain of DCLK1 altered its stability. Applying a combination of BioID-DCLK1 fusion proteins and phosphopeptide enrichment, Neurofilament medium polypeptide (NEFM), Filamin-A (FLNA), and Fatty acid synthase (FASN), mTOR, and CAMKK1 were found as possible DCLK1 interaction partners.

Novel Aspect

First phosphoproteomic analysis of neural progenitor cell differentiation identifying differential phosphorylation of DCLK1 and follow up by BioID fusion proteins.

References

1. Murry C., Keller G., *Cell*, 132, 661-680 (2008)
2. Miyamoto Y., Yamauchi J., Chan J.R., Okada A., Tomooka Y., Hisanaga S., Tanoue A., *J. Cell Sci.*, 120, 4355–4366 (2007)
3. Guardiola-Diaz H.M., Ishii A., Bansal R., *Glia*. 60, 476-486
4. Roux K.J., Kim D.I., Raida M., Burke B., *J. Cell Biol.* 196, 801-810 (2012)

TP-133 / ANALYSIS OF PROTEIN-PHENOLIC COMPOUND MODIFICATIONS USING ELECTROCHEMISTRY COUPLED TO MASS SPECTROMETRY

Constanze Kallinich - Simone Schefer - Sascha Rohn

Hamburg School of Food Science, University of Hamburg, Department Food Chemistry, Hamburg

Keywords: Electrochemical oxidation, phenolic compounds, adduct formation, EC-MS, phenol-protein interaction

Introduction

Posttranslational modifications of proteins are associated with influencing the essential structural and regulatory properties of each protein. For example, one of these post-translational modifications is the covalent addition of other food compounds to diverse side-groups of (food) proteins or peptides [1]. Influences on the biological properties such as the enzyme activity [2], the nutritional protein quality [3], and changes in the allergenicity of the proteins/protein derivatives, seem to be more or less obvious.

Methods

Electrochemical oxidation of phenolic acids (*e.g.*, chlorogenic acid) was performed using a preparative electrochemical thin-layer cell consisting of a boron-doped diamond working electrode. A solution of the phenolic acid was injected into the electrochemical cell and a constant potential was applied [4,5]. Detection of oxidation products was performed using an ESI-MS ion trap mass analyzer in negative ion mode. Subsequently, a second flow system was used to combine both analyte solutions, phenolic acid and mixture of tryptic peptides of a dietary protein *e.g.*, alpha-lactalbumin, to generate and detect potential adducts.

Results

After optimizing the reaction parameters for the four phenolic compounds selected, various oxidation products were detected and identified by means of ESI-MS. By connecting a flow system with a tryptic peptide solution of digested dietary proteins, oxidation products of the phenolic substances, and further reaction products of phenolic oxidation products and tryptic peptides were identified. It has been shown that an electrochemically induced oxidation of the model phenols is possible as well as a reaction of the oxidation products with selected, previously tryptically digested dietary proteins. A preparation of the reaction products in higher concentrations could be generated by preparative LC. Structural elucidation of the resulting products by means of NMR are carried out.

Conclusions

Electrochemistry coupled with mass spectrometry was successfully used to investigate the redox potential of dietary phenolic acids such as chlorogenic acid. It was possible to create optimized conditions for the oxidation of chlorogenic acid to generate a wide variety of oxidation products. Electrochemically oxidation of chlorogenic acid, caffeic acid, ferulic acid, and sinapinic acid generated a variety of reactive intermediates. Feasibility of this technique for the research topic was shown.

Novel Aspect

The young field of electrochemical coupled mass spectrometry (EC-MS) is a suitable system for generating and identifying potential metabolites. Target compounds can be identified directly by mass spectrometry. After separation of generated metabolites, structural elucidation can occur in following steps.

References

1. Le Bourvellec, C.; Renard, C., *Crit. Rev. Food Sci. Nutr.*, 52, 213–248 (2012)
2. Rohn, S.; Rawel, H. M.; Kroll, J., *J. Agric. Food Chem.*, 50, 3566–3571 (2002)
3. Petzke, K. J.; Schuppe, S.; Rohn, S.; Rawel, H. M.; Kroll, J., *J. Agric. Food Chem.*, 53, 3714-3720 (2005)
4. Büter, L.; Frensemeier, L. M. ; Vogel, M. ; Karst, U. *J. Chromatogr. A*, 1479, 153-160 (2017)
5. Jahn, S.; Faber, H. ; Zazzeroni, R. ; Karst, U., *Rapid Commun. Mass Spectrom.*, 26, 1415-1425 (2012)

TP-134 / A NEW PHOSPHOPROTEOME ENRICHMENT STRATEGY IN SMALL-SCALE USING PHOSPHO-SPECIFIC ANTIBODIES WITH ONLINE MHFER-NANOLC-ESI-MS/MS

Sun Young Lee - Kwonseong Kim - Dukjin Kang

Korea Research Institute of Standards and Science, Division of Chemical and Medical Metrology, Daejeon

Keywords: Phosphoprotein enrichment, phospho-specific antibodies, mHFER, nLC-ESI-FT orbitrap-MS/MS

Introduction

Phosphoprotein is involved in cellular processes such as proliferation, cell signaling, and apoptosis. [1-3] Mass spectrometry is commonly used to identify both phosphorylated protein(s) and its site(s). However, phosphoproteins in MS analysis have several limitations due to low abundant and low ionization efficiency caused by negative charge of phosphates in a positive mode. [4-5] Herein, we introduce a MW sorting-based new strategy for phosphoproteomics.

Methods

The tryptic peptides from MCF7 were mixed with phospho-antibodies (pAbs, above 150 kDa), and followed by which the mixtures introduced into the mHFER (50-100 kDa in MW-cutoff) with continuous flow rate of 5 μ L/min. During breakthrough run, ordinary peptides having no affinity with pAb are first washed out through membrane wall of HF. The remaining peptides bound to pAb can be simultaneously exited via online tryptic digestion in mHFER for next nLC-ESI-MS/MS run.

Results

For mHFER-based phosphopeptide enrichment, we optimized the experimental conditions according to the kind of pAbs, reaction time between peptides and pAbs, reaction temperature, reaction ratio, and the length of mHFER. From our experimental results, the enrichment yield of phosphopeptides was dramatically enhanced when used the mixture of phospho-serine, -threonine and -tyrosine antibodies at the ratio of 1:1 (proteins/pAbs, w/w) reacted for 1 hours at 4°C. Finally, we found total 11282 phosphorylation site from 5673 phosphopeptides of MCF7 cells utilizing online mHFER-2D-nLC-ESI-MS/MS and confirmed that the pAbs-based enrichment coupled with online mHFER can be the best way with the small amount of proteins (~10 μ g) for phosphoproteomics than conventional methods [e.g., immobilized

metal affinity chromatography (IMAC), titanium dioxide (TiO₂), and filter aided sample preparation (FASP)] with a large amount (above 500 ug).

Conclusions

We found that the developed pAbs-based online mHFER is highly selective, precise, and suitable for phosphopeptide enrichment with the small amount of proteins, compared to that of conventional methods (IMAC, TiO₂, and FASP). Additionally, our phosphopeptide enrichment strategy with online mHFER is a promising tool for phosphoproteomics, and for the development of therapeutic agents.

Novel Aspect

The developed a mHFER-based phosphoproteome enrichment platform is one of advanced strategy for understanding the biological functions of protein PTMs study.

References

1. Imamura H., Sugiyama M., Wakabayashi M., Ishihama Y., J. Proteome Res. 13, 3410 (2014).
2. Boehm M. E., Adlung L., Schilling M., Roth S., Klingmüller U., J. Proteome Res. 13, 5685 (2014).
3. Glibert P., Meert P., Steendam V., Nieuwerburgh F. V., Coninck D. D., Martens L., Dheanens M., Deforce D., J. Proteome Res. 14, 839 (2015).
4. Kweon H. K., Håkansson K., Anal. Chem. 78,1743 (2006).
5. WangX., Stewart P. A., Cao Q., Sna Q. -X. A., ChungL. W. K., Emmett M. R., Marshall A. G., J. Proteome Res. 10, 3920 (2011).

TP-135 / HISTONE MIDDLE-DOWN ANALYSIS USING BY MULTIPLE DISSOCIATION METHODS

Takeshi Kawamura⁽¹⁾ - ***Hiroki Ueda***⁽²⁾ - ***Daisuke Higo***⁽³⁾ - ***Yoko Chikaoka***⁽¹⁾ - ***Kazuki Yamamoto***⁽¹⁾ - ***Lee Sujin***⁽¹⁾

Isotope Science Center, The University of Tokyo, Tokyo (1) - Research Center for Advanced Science and Technology, The University of Tokyo, Tokyo (2) - Thermo Fisher Scientific, Tokyo (3)

Keywords: Histone modifications, epigenome, ETD, UVPD

Introduction

Mass spectrometry is a powerful tool for protein post-translational modifications (PTMs) analysis. A state-of-the-art mass spectrometer distinguishes very small mass differences. Histone has many PTM sites and their theoretical PTM combinations can be above trillion. Semi-bottom-up and middle-down proteomics are used for PTMs identification or decision of their combinations. Here, we used middle-down proteomics and compared various fragmentation methods.

Methods

Purified histone or histone H3 HPLC fraction were digested with AspN to obtain 49 or 23 aa length of H3 or H4 tails. They were fragmented with CID, ETD, or UVPD method and were analyzed by Mascot database search for PTM identification. The identification number and the identified kinds of PTMs were compared. The search was performed after deconvolution into the masses of singly charged ions.

Results

Although all fragmentation methods enabled PTM identifications, ETD showed relatively accurate identification in the middle-down proteomics. However, it was difficult to differentiate citrullination from z, z+1, z+2 ions generated by hydrogen rearrangement in ETD. In addition, the molecules of molecular weights above 5 kDa with multiple charged ions showed deformed isotope pattern of distributions, which make it difficult to determine monoisotopic ions. We have developed an algorithm to determine

monoisotopic ions by matching to theoretical isotope distributions and compared with the Advanced Precursor Detection (APD) in the Orbitrap Fusion Lumos. Moreover, CID and UVPD were compared and the CID was suitable for shorter peptide analysis. Longer peptide analysis was better in UVPD and the UVPD generated relatively frequent b, y series ions.

Conclusions

Appropriate combination of ETD, CID, and UVPD and selection of the optimized fragmentation method according to molecular size and charge might improve the accuracy of histone PTMs analysis. Further, the PTM site locations are examined.

Novel Aspect

We have compared the fragmentation methods to reveal complicated histone PTMs.

References

1. Sidoli et al., Expert Review of Proteomics, 14, 617–626 (2017)
2. Yamamoto et al., Mass Spectrometry, 4, A0039 (2015)
3. Greer et al., Proteome Res., 17, 1138–1145 (2018)

TP-136 / INFLUENCE OF DIGESTION CONDITIONS ON POST-TRANSLATIONAL MODIFICATIONS IN MONOCLONAL ANTIBODIES

Maria Khalikova⁽¹⁾ - ***Alexia Ortiz***⁽²⁾ - ***Koen Sandra***⁽²⁾

Charles University, Faculty of Pharmacy, Hradec Kralové (1) - Research Institute for Chromatography (R. I. C.), President Kennedypark 26, Kortrijk (2)

Keywords: Monoclonal antibodies, post-translational modifications, trypsin digestion, low pH digestion, peptide mapping

Introduction

Post-translational modifications (PTMs) on monoclonal antibodies (mAbs) are potential critical quality attributes that define therapeutic efficacy and safety[1]. Mass spectrometry based peptide mapping (PM) is the gold standard methodology for the in-depth study of PTMs on mAbs. Artificial PTMs, however, may be induced during sample preparation for PM, *e.g.* asparagine deamidation [2]. The current work describes the effect of digestion pH and incubation time on Trastuzumab PTMs. Trastuzumab, marketed as Herceptin, is a humanized IgG1 mAb used in the treatment of HER2 positive breast cancer.

Methods

The regular trypsin digestion protocol for peptide mapping involved the use of RapiGest, dithiothreitol (DTT), and iodoacetamide (IAM) in Tris-HCl buffer at pH 7.5. The low pH protocol, based on a commercial kit from Promega, made use of guanidine-HCl, Tris(2-carboxyethyl)phosphine (TCEP), IAM, trypsin, and low pH resistant rLys-C. Peptide mapping was performed on an Agilent Technologies 1290 Infinity UHPLC and Agilent Technologies 6540 Accurate-Mass Quadrupole Time-of-Flight LC/MS. Peptides were separated on a Waters Acquity UPLC CSH 130 C18 (2.1 x 150 mm, 1.7 μm) analytical column operated at 60°C at a flow rate of 0,4 mL/min using 0,1% formic acid in water or acetonitrile.

Results

Post-translational modifications of Trastuzumab peptides produced under regular alkaline, and low pH digestion conditions were identified and compared. Aspartate isomerization, methionine oxidation, *N*-glycosylation, *N*-terminal cyclization, and asparagine deamidation were PTMs of primary interest. Digests with incubation periods of 2, 4, 6, and 16 hours were obtained and analyzed by LC-MS.

Results show the progression of artificial modifications, especially asparagine deamidation, using regular tryptic digestion conditions with longer incubation periods up to 16 hours at 37°C. On the contrary, the low pH digestion protocol showed less identified artificial modifications while providing longer peptides and higher sequence coverage. Moreover, modifications not associated with sample preparation were identified, such as oxidation, deamidation, glycosylation, and cyclization of certain peptides.

Conclusions

The accurate determination of PTMs is important to ensure the quality of therapeutic mAbs. Artificial IPTMs introduced by sample preparation using alkaline tryptic digestion can be reduced with shorter incubation times while providing sufficient peptide profile and protein sequence coverage. Acidic digestion effectively prevents modifications caused by sample preparation.

Novel Aspect

Post-translational modifications of Trastuzumab under different digestion conditions were identified and compared. Shorter incubation times for regular trypsin digestion or the use of low pH digestion conditions is advised.

References

1. Jenkins N., Murphy L., Tyther R. *Molecular Biotechnology* 39, 113-118 (2008).
2. Hao P., Ren Y., Datta A., Tam J. P., Sze S. K. *Journal of Proteome Research* 14, 1308-1314 (2015)

This work was supported by the STARSS project (Reg. No. CZ.02.1.01/0.0/0.0/15_003/0000465) co-funded by ERDF, SVV 260 412.

TP-137 / IDENTIFICATION AND CHARACTERIZATION OF O-XYLOSYLATION AND O-GLYCOSYLATION OF RECOMBINANT PROTEIN THERAPEUTICS BY MASS SPECTROMETRY

Vincent Larraillet

Roche Pharmaceutical Research and Early Development, Biochemical and Analytical Research, Roche Innovation Center Munich

Keywords: antibody fusion proteins, post-translational modification

Introduction

To connect the different functionalities of antibody fusion proteins in a single molecule G4S linker are often used. O-Xylosylation is a recently described post-translational modification (PTM) of a heterogeneous glycan structure on the serine of the G4S linker [1, 2]. Content, structure vary between different therapeutic proteins and expression hosts. O-Glycosylation is a known PTM with well described glycan structures [3, 4] and located on serine or threonine.

Methods

Intact measurement: Protein samples were N-deglycosylated and desalted by size exclusion chromatography. ESI mass spectra were recorded on a Q-TOF instrument (maXis, Bruker) equipped with a nano ESI source (TriVersa NanoMate, Advion).

Peptide Mapping: Protein samples were denatured and reduced, carboxymethylated and digested by trypsin. LC/ESI-MS/MS was done on an Orbitrap Fusion system (Thermo). Peptides were separated using RP HPLC (Waters NanoAcquity).

Results

Here we describe mass spectrometry based methods to detect, differentiate and locate O-Xylosylation and

O-Glycosylation in recombinant protein samples. We will show that intact measurement of deglycosylated antibody derivatives using ESI-TOF-MS is an easy and reliable method to quickly get an initial impression of present protein modifications which exist with significant amount. A major disadvantage of this method is that detected O-Xylosylation signals are often overlaid by other modification or adduct signals or detected peaks are not fully resolved. Measurements after deglycosylation and reduction can provide remedy regarding peak resolution issues. But, as we will show, LC/MS-MS peptide mapping using an Orbitrap Fusion is the method of choice to get a more detailed view on the modification, to locate the particular modification site in protein sequence and to make an exact estimation of the O-Xylosylation and O-Glycosylation content.

Conclusions

Detection and quantification of O-Xylosylation and O-Glycosylation in antibody fusion proteins using the combination of intact measurement and peptide mapping.

Novel Aspect

Detection of O-Xylosylation and O-Glycosylation of recombinant protein therapeutics.

References

1. Spencer D., et al., *Journal of Pharmaceutical Sciences*, 102, 3920-3924 (2013).
2. Wen D., et al., *Analytical Chemistry*, 85, 4805-4812 (2013).
3. Peter-Katalinić J., *methods in Enzymology*, 405, 139-171 (2005).
4. Wuhler M., et al., *Molecular & Cellular Proteomics*, 14, 1373-1384 (2015).

TP-138 / SPECIFIC MASS TAG OF METHYLATED AMINO ACID VALIDATES THE SUMOYLATION SITES BY HCD ON ORBITRAP MASS SPECTROMETER

Fu-An Li - Yu-Hsiang Cheng

Institute of Biomedical Sciences, Academia Sinica, Taipei

Keywords: SUMOylation, methylated, HCD, Orbitrap

Introduction

The identification of SUMOylation sites by mass spectrometry is not easy because of the large SUMO chain remnant attached to peptides after tryptic digestion. Therefore, the widely adopted method is still the site-directed mutagenesis. Here, we will demonstrate the utility of specific mass tag from the methylated terminal to facilitate and validate the identification of SUMOylation sites by HCD on Orbitrap mass spectrometer without mutagenesis.

Methods

The digested peptides are reductively methylated using formaldehyde-d₂ and subsequently subjected to LC/MS/MS analysis with LTQ-Orbitrap Elite. The digested peptides were analyzed in data-dependent mode by LC/MS/MS. The data-dependent acquisition method consisted of ten sets of HCD/CID scans. Here, the HCD spectrum is used for monitoring specific a1 ions, and the CID spectrum is used for peptide sequencing.

Results

The signal enhancement for a1 ions of dimethyl labeling peptides in the CID of Q-TOF instrument had observed and applied for identifying the N-terminal amino acid. (1-3) In this study, we identified the spectra of SUMOylated peptides by monitoring two a1 ions of HCD spectrum, one is from the SUMO chain remnant and the other is from the SUMOylated peptide. The optimal collision energy of HCD for monitoring two a1 ions simultaneously was 40%.

For peptides sequencing, the identification of the MS/MS spectra generated by the SUMOylated peptide using the Mascot database search is still difficult. To address this problem, the candidate spectra selected by specific mass tag were submitted to do the sequence match by the software "ChopNSpice". (4) The results showed the specific mass tag of methylated amino acid could dramatically narrow down the candidates of possible spectra, thus the number of identified spectra were more than the one without spectra selection before ChopNSpice search.

Conclusions

We use a di-SUMO-2 protein as a model to demonstrate the method and successfully applied to identify the SUMOylated sites on the p53 protein. The results show the methylated mass tag (two specific a1 ions) is suitable for screening SUMOylated spectra in complex protein mixture. This method is more sensitive than direct database search, and it also greatly reduced the false identification.

Novel Aspect

The method could dramatically narrow down the candidates of possible spectra and is more sensitive than direct database search for SUMOylation sites identification.

References

1. Hsu J. L., Huang S. Y., Shiea J. T., Huang W. Y., Chen S. H., J Proteome Res., 4, 101-8 (2005).
2. Huang S. Y., Hsieh Y. T., Chen C. H., Chen C. C., Sung W. C., Chou M. Y., Chen S. F., Anal Chem., 84, 4900-6 (2012).
3. Chicooree N., Griffiths J. R., Connolly Y., Smith D. L., Rapid Commun Mass Spectrom., 27, 2108-14 (2013).
4. Hsiao H. H., Meulmeester E., Frank B. T., Melchior F., Urlaub H., Mol Cell Proteomics., 8, 2664-75 (2009).

TP-139 / HOW DOES CE-MS COMPARE WITH LC-MS IN PTM ANALYSIS?

Klaus Faserl⁽¹⁾ - **Bettina Sarg**⁽¹⁾ - **Herbert Lindner**⁽¹⁾ - **Jim Thorn**⁽²⁾ - **Christopher Lößner**⁽³⁾ - **Stephen Lock**⁽²⁾

Medical University of Innsbruck, Division of Clinical Biochemistry, Innsbruck (1) - SCIEX, CE and Biopharma, Warrington (2) - SCIEX, CE and Biopharma, Darmstadt (3)

Keywords: PTM, CE-MS, phosphorylation, deamidation, isomerisation

Introduction

Post translational modifications (PTMs) are important indicators of change in cells. Two of the top four most abundant PTMs are deamidation and phosphorylation and the location of the phosphorylation site is an important factor in understanding the effect of this modification on the activity of proteins.

Results

Traditional LC-MS methods can sometimes struggle with the identification of positional isomers of phosphorylated peptides as they are identical in mass and may not separate by traditional reverse phase LC. Phosphorylated peptides are often polar, and elute early in reverse phase chromatography (especially when more than one phosphorylation site is present on the peptide) which can make them difficult to detect by trap-elute LC-MS. LC-MS methods also struggle to identify and quantify Aspartate and iso-Aspartate isomers (associated with deamidation) as they have the same mass and similar fragmentation patterns which often result in false positive identifications.

Conclusions

As CE separates analytes based on their charge and size (a different mechanism to LC) it can often overcome some of the separation challenges seen by LC. In this poster we will show how CESI-MS has been used to tackle these challenges and how it compares to LC-MS and find out how it has been applied to real

biological samples in proteomics research.

Novel Aspect

Separation of isobaric and/or co-eluting post-translational modifications by CE-MS

TP-140 / CHARACTERIZATION OF NIST MONOCLONAL ANTIBODY ON INTACT, SUBUNIT, AND PEPTIDE LEVEL WITH MONITORING OF CQA'S ON THE 6600

Chris Nortcliffe⁽¹⁾ - *Kerstin Pohl*⁽²⁾ - *Jörg Dojahn*⁽³⁾ - *Ferran Sanchez*⁽⁴⁾

SCIEX, Phoenix House, Lakeside Drive, Warrington (1) - SCIEX, Applications Support, Darmstadt (2) - Sciex, Market Manager, Darmstadt (3) - Sciex, Market Manager, Warrington (4)

Keywords: Mass Spectrometry, Biologics, Antibodies, MAM, Biopharmaceuticals

Introduction

Monoclonal antibodies (mAbs) can be incredibly complicated with hundreds of residues, varying degrees of glycosylation and PTMs which may be critical to quality. MS can uncover a wealth of information of these mAbs with a synergistic workflow used to maximise data quality and fidelity. Here BioPharmaView™ and the 6600 characterize the NIST mAb on an intact, subunit, and peptide level. Additionally the MAM workflow in BioPharmaView™ has been used to control CQAs.

Methods

All measurements were carried out in replicates on a 6600 Q-ToF equipped with IonDrive™ source and the high flow ExionLC™ system. For all analyses solvent A consisted water containing 0.1 % formic acid and Solvent B consisted acetonitrile containing 0.1 % formic acid. Intact and subunit analysis were performed with a C4 column with a flow rate of 200 µl/min at temperatures of 80°C and 60°C respectively. Peptide mapping was performed with a C18 column at 250 µl/min at 50°C.

Results

As complexity of biopharmaceuticals increases the need for accurate reliable data acquisition becomes ever more important. The NIST mAb has five major intact glycoforms G0F/G0F, G0F/G1F, G1F/G1F, G1F/G2F, and G2F/G2F representing one glycan on each heavy chain Fc region. The ratios of these various glycoforms are highly regulated and have been measured accurately after reconstruction with BioPharmaView™. Due to the number of overlapping features at high charges states features can be difficult to discern at the intact level, these features can become clearer at the subunit level. The glycoforms on the subunit level are clearer to distinguish due to the reduced complexity of a single glycan per fragment. The subunit fragments were matched with the BioPharmaView™ software providing the same glycan ratios as on the intact level. The NIST mAb was digested and analysed with the 6600 with a SWATH analysis. Data was analysed with BioPharmaView™ software with 100 % of the protein sequence observed at the MS level and validated via MS/MS fragmentation.

Conclusions

As biologics become increasingly heterogeneous the need for complete characterization across all points in upstream and downstream production increases. Sciex provides a complete workflow solution for LC-MS analysis of complex protein products with Multiple Attribute Monitoring of Critical Quality Attributes.

Novel Aspect

Complete characterization of NIST antibody and glycoforms with MAM analysis in BioPharmaView™ software.

References

1. Formolo T., Ly M., Levy M., Kilpatrick L., Lute S., Phinney K., Marzilli L., Brorson K., Boyne M., Davis D., Schiel J., State-of-the-Art and Emerging Technologies for Therapeutic Monoclonal Antibody Characterization Volume 2. Biopharmaceutical Characterization: The NISTmAb Case Study, 1-62 (2015)

TP-141 / THE PASEF METHOD ON A TIMS-QTOF MASS SPECTROMETER FOR HIGH SENSITIVITY PHOSPHOPROTEOMICS

Heiner Koch⁽¹⁾ - *Thomas Kosinsky*⁽¹⁾ - *Matt Willets*⁽²⁾ - *Markus Lubeck*⁽¹⁾ - *Scarlet Koch*⁽¹⁾ - *Oliver Raether*⁽³⁾ - *Gary Kruppa*⁽⁴⁾ - *Anjali Alving*⁽⁵⁾

Bruker Daltonik GmbH, Proteomics, Bremen (1) - Bruker Daltonics inc, Application, Billerica (2) - Bruker Daltonik GmbH, Development, Bremen (3) - Bruker Daltonics Inc, Proteomics, Billerica (4) - Bruker Daltonics Inc, Billerica, USA (5)

Keywords: TIMS, PASEF, Phosphopeptides, isobaric

Introduction

Analysis of protein phosphorylation by mass spectrometry serve as a powerful tool for biologists and clinical researchers. However, the abundance differences as well as the variety of phosphorylation sites are pushing the need for instruments with higher sensitivity and selectivity. Here we evaluate the performance of a dual trapped ion mobility separation (TIMS) - QTOF that provides an additional dimension of separation and higher peak capacity.

Methods

Phosphopeptides were enriched from 50-200 µg of proteolytic digests, separated using 15 to 90 min nano-HPLC gradients (nanoElute, Bruker Daltonics) and analyzed on a timsTOF Pro instrument using the Parallel Accumulation Serial Fragmentation (PASEF) method (Bruker Daltonics). Post-processing analyses were performed in PEAKS (Bioinformatics Solutions Inc.) and MaxQuant (Jürgen Cox, Max Planck Institute of Biochemistry).

Results

Applying the PASEF method on the timsTOF Pro instrument enables on average the selection of 12 precursors within 100 ms. In addition to the high speed, sensitivity in the timsTOF Pro is achieved by a close to 100 % duty cycle of the dual TIMS cell design, cleanup of MS/MS spectra by ion mobility and time and space focusing in the TIMS analyzer. Using the combination of speed and sensitivity in PASEF method we identified more than 66,000 peptide spectrum matches (PSMs) resulting in around 25,000 unique peptides and 16,000 unique phosphopeptides within a 90 min gradient run. In addition to m/z separation, phosphopeptides are separated by space to charge, which increases the peak capacity with the most pronounced effect on very short gradients. Together, the high speed of PASEF phosphopeptide separation in the m/z to Ω/z space, combined lead to 7,000 and more than 10,000 unique phosphopeptides identified at 15 min and 30 min gradient times respectively, at a speed of up to 500 unique phosphopeptides/min. Some of those separated phosphopeptides were isobaric.

Conclusions

TIMS-PASEF allows to increase phosphopeptide analysis depth and selectivity.

Novel Aspect

High speed and sensitivity provided by PASEF for fast and comprehensive phosphoproteome analysis.

TP-142 / THE CITRULLINE EFFECT IN THE DISSOCIATION OF DEIMINATED PEPTIDES

Arnold Steckel⁽¹⁾ - *Katalin Uray*⁽²⁾ - *Dávid Papp*⁽³⁾ - *Gitta Schlosser*⁽⁴⁾

Doctoral School of Chemistry, ELTE Eötvös Loránd University, Department of Organic Chemistry, Budapest (1) - MTA-ELTE Research Group of Peptide Chemistry, ELTE Eötvös Loránd University, Budapest (2) - Department of Organic Chemistry, ELTE Eötvös Loránd University, Budapest (3) - Department of Analytical Chemistry, ELTE Eötvös Loránd University, Budapest (4)

Keywords: citrulline, citrullination, citrulline effect, tandem mass spectrometry, peptide sequencing

Introduction

Protein deimination is a post-translational modification of proteins converting certain arginines to citrullines. Overcitrullination could be associated with severe pathological conditions. Mass spectrometric analysis of the modified proteins is hindered by several problems. Based on our tandem mass spectrometric experiments, we report here a new and pronounced amino acid effect.

Methods

Model peptides containing either one or two citrulline residues as well as their native, arginine-containing variants were synthesized. Their solutions were injected directly to the electrospray source of a high accuracy and resolution quadrupole-time of flight instrument and were analysed by tandem mass spectrometry using low-energy collision induced dissociation.

Results

A preferential cleavage site was observed at the amide bond of citrulline and any other following alpha amino acid (Cit-Aaa) yielding intensive complementary b and/or y type ions. If citrulline was located at the C-termini of the peptides, an enhanced cleavage of Aaa-Cit amide bond was perceived instead. These phenomena were described by us for the first time [1] and was coined citrulline effect, analogously to ornithine effect reported previously [2]. Loss of isocyanic acid from the citrulline side chains observed earlier [3] was also significant and effective for confirming the site of modification.

Conclusions

Fragmentation behaviour of citrulline containing peptides is similar to ones not containing them, however, significant intensity differences can be observed due to the newly identified citrulline effect. The loss of isocyanic acid is also a pronounced fragmentation pathway. Our results are in concordance with the mobile proton hypothesis.

Novel Aspect

The citrulline effect facilitates the MS/MS confirmation of citrulline residues in a peptide sequence, besides the neutral loss of isocyanic acid.

References

1. Steckel A., et al., Rapid Commun. Mass Spectrom., accepted article (2018)
2. McGee WM, McLuckey SA., J. Mass Spectrom., 48, 856-861 (2013)
3. Hao G. et al., J. Am. Soc. Mass Spectrom., 20, 723–727 (2009)

Acknowledgement

Gitta Schlosser acknowledges the support of the MTA Premium Post-Doctorate Research Program of the Hungarian Academy of Sciences (HAS, MTA).

TP-143 / CHEMICAL ACETYLATION OF HISTONE LYSINE RESIDUES AT THE PROTEIN LEVEL USING ACETIC ANHYDRIDE

Alienke van Pijkeren - *M. D. Kwiatkowski* - *Rainer Bischoff*

Keywords: histones, chemical acetylation, isotopically-labeled acetic anhydride, sodium borate

Introduction

Age-related diseases are associated with alterations in energy metabolism and with chronic inflammation. Research indicates that changes in energy metabolism are closely linked to gene expression via changes in histone acetylation profiles, possibly leading to a pro-inflammatory phenotype [1]. To investigate the histone acetylation dynamics by mass spectrometry, an efficient chemical acetylation reaction of unmodified histone lysine residues is necessary.

Methods

Histones derived from RAW264.7 and 16HBE cells were chemically acetylated using stable isotope labeled acetic anhydride in combination with different buffer systems; ammonium bicarbonate (ABC), triethylammonium bicarbonate (TEAB) or sodium borate (SB). After in-solution tryptic digestion, peptides were analyzed by liquid chromatography tandem-mass spectrometry (LC-MS/MS). Chemical acetylation efficiency was assessed by using Xcalibur software.

Results

LC-MS/MS analysis resulted in identification of histone proteins (H1, H2A, H2B, H3 and H4) with a high sequence coverage, between 86% and 98%. The average acetylation efficiency of 14 selected lysine residues of the core histone proteins H2A, H3 and H4 was significantly increased by using isotopically-labeled acetic anhydride in combination with either the TEAB and SB buffer (on average 99.6%) compared to the ABC buffer (94%). For some acetylation sites (for example H3K14, H3K56 and H4K20), no considerable differences in the acetylation efficiency were observed between the three different buffer systems. For other acetylation sites (for example H3K18, H3K79 and H4K20), the acetylation efficiency achieved by the TEAB and SB buffer was significantly increased compared to the ABC buffer. Complete acetylation of the histone 4 4-17 peptide, which contains four lysine residues, was significantly increased as well by using the TEAB and SB buffer (95%) compared to the ABC buffer (84%).

Conclusions

Chemical acetylation of unmodified histone lysine residues by using isotopically-labeled acetic anhydride in combination with TEAB or SB provided, in contrast to ABC, a sufficient acetylation efficiency (>95%). This acetylation method can be used to study the histone acetylation turnover in age related diseases. Furthermore, the method can be used to study the effects of histone deacetylase (HDAC) inhibitors on the acetylation status of specific lysine residues.

Novel Aspect

Acetic anhydride in combination with SB results in an almost complete chemical acetylation of histone lysine residues. SB outcompetes the more commonly used ABC buffer.

References

1. Peleg S., Feller C., Ladumer AG., Imhof A., Trends in Biochemical Sciences, 41, 700-711 (2016).

TP-144 / HIGH RESOLUTION MASS SPECTROMETRY BASED APPROACH FOR THE CHARACTERIZATION OF PLASMA PROTEINS OXIDATIVE MODIFICATIONS BY ROS

Alexandra Vasilyeva⁽¹⁾ - *Lyubov' Yurina*⁽¹⁾ - *Maria Indeykina*⁽²⁾ - *Anna Bugrova*⁽³⁾ - *Alexey Kononikhin*⁽²⁾ - *Mark Rosenfeld*⁽¹⁾ - *Evgeny Nikolaev*⁽⁴⁾

N. M. Emanuel Institute of Biochemical Physics, RAS, Laboratory of Thermodynamics of Biosystems, Moscow (1) - N. M. Emanuel Institute of Biochemical Physics, Institute for Energy Problems of Chemical

Physics, RAS; Moscow Institute of Physics and Technology, Laboratory of mass spectrometry of biomacromolecules, Moscow (2) - N. M. Emanuel Institute of Biochemical Physics, RAS, Laboratory of Neurochemistry, Moscow (3) - N. M. Emanuel Institute of Biochemical Physics, Institute for Energy Problems of Chemical Physics, RAS; Skolkovo Institute of Science and Technology, Laboratory of mass spectrometry of biomacromolecules, Moscow (4)

Keywords: Proteins, ROS, PTM, fibrinogen, fibrin stabilizing factor, plasminogen, HPLC-MS/MS

Introduction

Free radical oxidation in the human body causes the accumulation of oxidized forms of proteins, which have an impaired function and can also be used as the markers of many pathological conditions, induced diseases such as cancer, nephropathy, and cardiovascular diseases. High resolution mass spectrometry approach is promising for characterization of oxidative protein modifications by reactive oxygen species (ROS) [1].

Methods

Protein isolation, analyses of enzyme activity and the polymerization of fibrinogen were accomplished as described in [2,3,4]. After the oxidation of samples by ozone or HOCl, the proteins were hydrolyzed, using a suitable technique for each protein. HPLC-MS/MS experiments were performed on an Agilent 1100 nanoLC (Agilent Technologies, USA) coupled to 7T LTQ-FT Ultra (Thermo, Germany). PTM searching was carried out by PEAKS Studio software (v. 8.5, BSI, Canada).

Results

The site-specific oxidative modifications in the proteins (fibrinogen, plasminogen, fibrin-stabilizing factor) under induced oxidation in model systems with various oxidants were analyzed by HPLC-MS/MS method.

Oxidative modifications have been found in some functionally significant amino acid residues, as well as in sites whose function remains unclear and should be further investigated.

The variation of the coverage in native and oxidized samples is compared. These results obtained were also used to analyze the effect of oxidation on the structure of the protein molecules [5].

Interestingly, along with the sulfur-containing residues of Met and Cys, aromatic residues of Tyr and Trp, polar residues of Glu, Lys, and Asp and non-polar Pro residues, are also targets for the oxidants.

The information obtained by biochemical methods in complex with HPLC-MS/MS data, make it possible to suggest a mechanism for the functional property loss of proteins under oxidation.

Conclusions

The data obtained make it possible to evaluate the availability of amino acid residues of proteins for oxidants under conformational changes in enzyme molecules, to reveal oxidation sites and the degree of oxidation of amino acid residues by using different concentrations of oxidants. Moreover it is one may suggest the mechanism for the functional property loss of proteins under oxidation [6].

Novel Aspect

High resolution mass spectrometry based approach was developed for the characterization of plasma proteins oxidative modifications by ROS.

References

1. I. Y. Petrushanko, V. M. Lobachev, A. S. Kononikhin and etc., PLoS ONE (2016).
2. N. J. White, Yi Wang, Xiaoyun Fu and etc., Free Radic. Biol. Med. 96, 181-189 (2016)
3. L. Lorand, R. Credo, T. Janus, Methods Enzymol. 80, 333-41 (1981).
4. D. G. Deutsch, E. T. Mertz, Science (Washington, D. C.) 170, 1095-1096 (1970)
5. A. Vasilyeva, A. V. Bychkova, A. E. Bugrova, M. I. Indeykina and etc., Dokl. Biochem. Biophys. 472, 40-43 (2017).
6. M. A. Rosenfeld, A. D. Vasilyeva, L. V. Yurina and etc., Free Radic. Res. 52, 14-38 (2018).

TP-145 / MASS SPECTROMETRY-BASED INVESTIGATIONS OF COLD PHYSICAL PLASMA-INDUCED PTM'S IN PEPTIDES

Sebastian Wenske - Jan-Wilm Lackmann - Giuliana Bruno - Sander Bekeschus - Klaus-Dieter Weltmann - Thomas von Woedtke - Kristian Wende

INP Greifswald e. V., ZIK plasmatis, Greifswald

Keywords: Post-translational modification (PTM), Cold atmospheric plasma (CAP), Mass spectrometry, Peptides

Introduction

Oxidative PTM's of proteins and peptides play a role in many signaling pathways with potentially drastic changes in structure and activity of this molecule by distinct chemical groups [1]. Creating a library of such PTM's facilitates global mapping of PTM's in mammalian tissue and ultimately foster understanding of complex redox controlled processes, *e.g.* wound healing. Here, cold atmospheric plasma (CAP) was used as means to simulate oxidative conditions on various peptides.

Methods

Different CAP sources were used for a determined modification of the model peptides Bradykinin, Angiotensin, and the artificial peptide LYTF AHD. The resulting products were analyzed by LC-HRMS/MS-measurements (Q Exactive™) by direct infusion ESI or by LC-ESI coupling using a PepMap column. A workflow for the identification of the modified peptides with various MS software solutions was developed by using different algorithm for an optimal result validation.

Results

Using manual annotation, a time dependent increase in number and quantity of PTM's was observed for the plasma treatment. Single and double oxidation and hydroxylation were identified in high intensities. For Bradykinin, the observed adducts were compared with literature [2] showing that only up to five oxygen adducts were found after kINPen treatment instead of 10 as described for the COST-jet. In contrast, the synthetic peptide carried up to 10 oxidation adducts after 5 minutes of plasma treatment. The two amino acids phenylalanine and tyrosine were oxidized particularly efficient. Investigation of additional modification patterns are currently under way. After testing different evaluation softwares for a manual identification of the PTM's, the next step will be the optimization of a semi-automatic workflow for PTM annotation and quantification by using isotopic labeled peptides.

Conclusions

Cold physical plasma can act as a rich source for oxidative PTM's of a peptide, with numerous modifications detected by LC-ESI-MS in the model used. Plasma sources differ in modification efficacy, *e.g.* the He-COST-Jet is more effective than Ar-kINPen, due to the presence of atomic oxygen. Using cold plasma, an oxidative and nitrosative focused PTM-library for peptides and amino acids is under construction and should ultimately forecast PTMs in complex samples.

Novel Aspect

Usage of cold physical plasma as a tool to create an oxidation focused PTM library using relevant model peptides with the ultimate goal of predicting PTMs in complex samplers.

References

1. Karve T. M., Cheema A. K., J. Amino Acids, 2011, 1-13 (2011).
2. Verlackt C. C. W., Van Boxem W., Dewaele D., Lemièrre F., Sobott F., Benedikt, J., Neyts E. C., Bogaerts A., J. Phys. Chem., 121, 5787-5799 (2017).

TP-146 / NEXT-GENERATION PROTEOMICS INSIGHTS INTO BACILLUS CEREUS VIRULOME HIGHLIGHTS THE KEY ROLE OF METHIONINE OXIDATION OF EXOPROTEINS

Beatrice Alpha-Bazin⁽¹⁾ – ***Jean-Paul Madeira***⁽²⁾ - ***H. Omer***⁽²⁾ - ***Jean Armengaud***⁽¹⁾ – ***Catherine Duport***⁽¹⁾

Cea, Drf/Joliot/Dmts, Bagnols-Sur-Ceze, France (1) - UMR0408/SQPOV, Avignon Universite, INRA, Avignon, France

Keywords: bacillus cereus, shotgun proteomics, methionine oxidation, comparative exoproteome

Introduction

Bacillus cereus is a food-borne pathogen responsible of toxi-infections in humans. During the course of infection the survival and growth of the bacteria depend on its capacity to maintain an intracellular redox state in response to environmental changes. Its virulence depends on the production of virulence factors secreted in the extracellular milieu. In order to determine the key factors that contribute to redox homeostasis and virulence, we have performed high-throughput shotgun proteomics.

Methods & Results

In the studies cellular proteome and exoproteome time courses were monitored by nanoLC-MS/MS. When required comparisons were performed between wild-type and the mutant of interest cells. The roles of the redox balance sensor Rex [2] and the organic hydroperoxide stress resistance system OhrRA [1] in the machinery that maintains intracellular redox state was deciphered by next-generation proteomics. Their role in toxinogenesis was highlighted through a detailed analysis of exoproteome [1,2]. Exoproteome analysis also highlighted that toxins and other extracellular virulence factors contained oxidized methionine residues[3]. We showed that oxidation of Met residues in extracellular proteins is regulated by methionine sulfoxide reductase AB (MsrAB), a key factor of the antioxidant system, suggesting that secretion of virulence factors could contribute to maintain intracellular redox homeostasis [4]. Furthermore, interactions of bacteria genomic information and their recently acquired plasmids were reported once again through the proteomics viewpoint [4]. We also contributed to identify by proteogenomics a new key virulence factor that we named EntD which was revealed to be pivotal for pathogenesis[5], meanwhile it was missed by automatic annotator software.

Conclusions

Bacillus cereus is well able to survive to various environment conditions. The bacteria exploits its metabolic flexibility to grow and produce extracellular toxins. We have shown that Met residues in exoproteins and their recycling by MsrAB serve as an antioxidant strategy. This project illustrates the potential of next-generation proteomics in microbiology.

References

1. Clair G., Lorphelin A., Armengaud J., Duport C., *J. of Proteomics*94, 527-539 (2013).
2. Laouami S., Clair G., Armengaud J., Duport C., *PLoS One* 9, e107354 (2014)
3. Madeira JP., Alpha-Bazin BM., Armengaud J., Duport C., *Front Microbiol.* 6, 1-15 (2015)
4. Madeira JP., Alpha-Bazin BM., Armengaud J., Duport C., *Front Microbiol.* 8, 1-16 (2017)
5. Madeira JP., Omer H., Alpha-Bazin BM., Armengaud J., Duport C., *J Proteomics* 146, 25-33 (2016)
6. Omer H., Alpha-Bazin BM, Brunet JL., Armengaud J., Duport C., *Front Microbiol* 6, 1-15 (2015)

TP-147 / IDENTIFYING THE DISULFIDE BOND PATTERN AND ECM INTERACTION SITE OF TGFBI

Marie V. Lukassen⁽¹⁾ - ***Carsten Scavenius***⁽²⁾ - ***Ida B. Thøgersen***⁽²⁾ - ***Jan J. Enghild***⁽²⁾

Aarhus University, Interdisciplinary Nanoscience Center (iNANO), Aarhus, Denmark (1) - Aarhus University,

Keywords: Disulfide, Cysteinylation, RP-HPLC, ECM, Cross-link

Introduction

Transforming growth factor β -induced protein (TGFBIp) is an extracellular matrix (ECM) protein composed of an NH₂-terminal cysteine-rich domain (CRD) and four fasciclin-1 (FAS1-1 to FAS1-4) domains. Previous studies have shown that TGFBIp is cross-linked to type VI and type XII collagen through a reducible bond [1,2]. TGFBIp contains 11 cysteine residues and is thus able to form five disulfide bonds, leaving a single cysteine residue available for the cross-link.

Methods

Recombinant human TGFBIp was digested with Glu-C and trypsin to get cleavages between all cysteine residues except Cys84 and Cys85. The generated peptides were separated by strong cation exchange (SCX) followed by reversed-phase high-performance liquid chromatography (RP-HPLC) and analyzed by mass spectrometry (MS) and NH₂-terminal sequencing. For verification, the disulfide bond pattern of TGFBIp from human cornea was also analyzed by MS.

Results

The NH₂-terminal CRD contains six cysteine residues, and one of these (Cys65) was identified as the candidate for the reducible cross-link between TGFBIp and the ECM molecules type VI and type XII collagen. This residue was not engaged in any disulfide bond, but was modified by a cysteinylation. The cysteinylation can be exchanged with a cysteine residue of another protein forming an intermolecular cross-link. In addition, the CRD contained two intradomain disulfide bonds (Cys49-Cys85 and Cys84-Cys97) and one interdomain disulfide bond to the FAS1-2 domain (Cys74-Cys339). The cysteine residues in the FAS1-3 domain (Cys473 and Cys478) were shown to form an intradomain disulfide bond. Finally, an interdomain disulfide bond between the FAS1-1 and FAS1-2 domains (Cys214-Cys317) was identified. The disulfide bond pattern was verified in human corneal tissue and Cys65 was also found in a cysteinylated form.

Conclusions

The results contradict the predicted assignment of the CRD as an emilin domain [3] and we have renamed the domain to the cysteine-rich domain of periostin and TGFBIp (CROPT). The interaction to the ECM is through Cys65 of the CROPT domain. The interdomain disulfide bonds indicate that the NH₂-terminus of TGFBIp adopts a compact globular fold. This together with the disulfide bond pattern has recently been validated by a crystal structure of TGFBIp [4].

Novel Aspect

We identified Cys65 of the renamed NH₂-terminal CROPT domain to be the residue forming a reducible bond to other ECM molecules such as type VI and type XII collagen.

References

1. Runager K., Klintworth G. K., Karring H., Enghild J. J., *Biochemistry*, 43, 2821-2827 (2013)
2. Hanssen E., Reinboth B., Gibson M. A., *Journal of Biological Chemistry*, 278, 24334-24341 (2003)
3. Callebaut I., Mignotte V., Souchet M., Mornon J., *Biochemical and Biophysical Research Communications*, 300, 619-623 (2003)
4. Garcia-Castellanos R., Nielsen N. S., Runager K., Thøgersen I. B., Lukassen M. V., Poulsen E. T., Goulas T., Enghild J. J., Gomis-Ruth F. X., *Structure*, 25, 1740-1750 (2017)

TP-148 / INCREASING VARIANTS IN HYDROPHOBIC INTERACTION CHROMATOGRAPHY DURING STABILITY STUDY OF HUMANISED MONOCLONAL ANTIBODY IDENTIFIED BY LC-MS ESI Q-TOF TOP DOWN, MIDDLE DOWN AND BOTTOM UP ANALYSES

Martin Pauers - Katharina List - Dorothee Wetzel

Boehringer Ingelheim Pharma GmbH & Co KG, Pharmaceutical Development Biologicals, Biberach an der Riss, Germany

Keywords: Biotherapeutics, humanized monoclonal antibodies, stability studies, hydrophobic interaction chromatography

Introduction

Development of therapeutic mAbs includes real time and accelerated stability studies to address critical Quality Attributes. Typically, these include *e.g.* protein structure, formation of particles as well chemical modifications of amino acid residues. This case study presents results from LC-MS analyses of samples from a stability study of a humanised IgG, which showed a storage temperature/ duration dependent increase of unknown variants within HIC analyses.

Methods

Samples stored at -70°C/ +5°C/ +25°C for 12 months and +30°C for 6 months were analysed by LC-MS ESI Q-TOF using a Waters UPLC I-Class coupled to a Waters Synapt G2-S. Performed analyses covered top down and multiple middle down analyses of IgG and its subunits using reducing conditions and/or IdeS treatment. Further analyses included bottom up analyses like tryptic and chymotryptic peptide mapping and determination of free thiol levels and disulfide mapping.

Results

LC-MS top down analyses resulted in detection of slightly increased protein fragment levels for samples stored at +25°C and +30°C, as well as detection of an additional hydrophilic light chain species without mass shift in middle down analyses. Bottom up analyses of these samples (+25°C/ +30°C) by peptide mapping resulted in slight increase of methionine oxidation, asparagine deamidation and aspartate isomerisation of DG motifs. Free thiol levels as well as scrambled disulfides in disulfide mapping showed slight increases at +25°C and +30°C. In-depth data evaluation of peptide mapping data indicated high levels of three additional retention time-shifted peptide species. These peaks were further investigated by LC-MS including electron transfer dissociation resulting in the identification of two discrete non-DG aspartate isomerisation sites in the variable domain of the light chain and one isomerisation site in the hinge region of the heavy chain.

Conclusions

LC-MSETD analyses of samples stored at +25°C/ +30°C detected elevated isomerisation levels of three non-DG motif aspartates. Two of them are located in the CDR1 and CDR2, therefore impact on the target binding was suggested. Analyses like DLS, fluorescence spectroscopy, RALS, ITF, IR and UV absorbance detected no major differences. Thus only in-depth LC-MS analyses revealed the root cause of additional HIC peak. Identified isomerisation sites were assessed as cQAs.

Novel Aspect

Unexpected high isomerisation levels for three Aspartates during the accelerated stress study have been identified.

TP-149 / HIGH RESOLUTION MASS SPECTROMETRY CHARACTERIZATION OF THE OXIDATION PATTERN OF METHIONINE AND CYSTEINE RESIDUES IN HUMAN MITOCHONDRIA VOLTAGE-DEPENDENT ANION SELECTIVE CHANNEL ISOFORMS

Rosaria Saletti⁽¹⁾ - *Maria G. G. Pittalà*⁽²⁾ - *Pierpaolo Risiglione*⁽²⁾ - *Vincenzo Cunsolo*⁽¹⁾ - *Vito De Pinto*⁽²⁾ - *Salvatore Foti*⁽¹⁾

Department of Chemical Sciences, Organic Mass Spectrometry Laboratory, University of Catania, Viale A. Doria 6, 95125 Catania, Italy (1) - Department of Biomedicine and Biotechnology, Section of Biology and Genetics, Viale A. Doria 6, 95125 Catania, Italy (2)

Keywords: Voltage dependent anion selective channel isoforms; Cysteine over-oxidation; Mitochondrial intermembrane space; High Resolution mass spectrometry

Introduction

Voltage-dependent anion selective channels are a quantitatively relevant component of the mitochondrial proteome. They are located in the outer membrane, which separates the cytosol from the inside of the organelle [1] and are responsible for the exchange of ions and small molecules in and out of mitochondria. Since the mitochondrial intermembrane space has an oxidative potential, we have investigated the oxidation state of cysteine and methionine residues of hVDACs.

Methods

Intact mitochondria, extracted from HAP1 cell culture, were reduced with DTT, alkylated by addition of iodoacetamide and then lysed. After centrifugation, the supernatant containing the mitochondrial proteins was loaded onto a hydroxyapatite column. The hVDACs enriched eluate obtained was desalted, in-solution digested using porcine trypsin or chymotrypsin, and the resulting peptide mixtures analyzed by UHPLC/High Resolution nanoESI-MS/MS.

Results

Our results demonstrate that all the three mitochondrial hVDACs, in physiological state, contain all the methionines partially oxidized to methionine sulfoxide. Furthermore, in the isoform 1 Cys127 is detected exclusively in the oxidized form of sulfonic acid, while the Cys232 is observed totally carboxyamidomethylated. In the hVDAC2, the cysteine residues 8, 13, 47, 133, 210, and 227 are only carboxyamidomethylated, Cys138 was identified trioxidized, and the cysteines 76 and 103 are partially in the form of sulfonic acid. Finally, also for the isoform 3 the results confirmed the presence of some cysteines (2, 8, and 229) reduced and the partial oxidation of Cys36 and Cys65. The peculiar behavior of Met and Cys residues of hVDACs may be related with the accessibility of the protein to a strongly oxidizing environment and may be connected with the regulation of the activity of this trans-membrane pore proteins.

Conclusions

In this work, adopting the same procedure developed for the structural investigations of rVDACs [2,3], a detailed profile of the oxidation state of methionine and cysteine residues have been obtained in all the three human isoforms. Our results demonstrate that the mitochondrial hVDACs, in physiological state, contains methionines oxidized to methionine sulfoxide and cysteines in the form of sulfonic acid.

Novel Aspect

The structural features elucidated by the present work may be helpful for a better understanding of the functional role of these proteins in the cell.

References

1. Benz R., *Biochimica Biophysica Acta*, 1197, 167-196 (1994).

2. Saletti R, Reina S, Pittalà MGG, Belfiore R, Cunsolo V, Messina A, De Pinto V, Foti S., *Biochimica Biophysica Acta*, 1859, 301-311 (2016).
3. Saletti R., Reina S., Pittalà M. G. G., Magrì A, Cunsolo V., Foti S., De Pinto V, *Biochimica Biophysica Acta-Bioenergetics*, (2018) doi.org/10.1016/j.bbabi.2018.06.007

TP-150 / UNDERSTANDING EPIGENETIC CONTROL OF HISTONE MODIFICATIONS DURING DROUGHT ACCLIMATION IN SORGHUM USING TOP-DOWN MASS SPECTROMETRY

Mowei Zhou⁽¹⁾ - **Neha Malhan**⁽¹⁾ - **Amirhossein Ahkami**⁽¹⁾ - **Kristin Engbrecht**⁽¹⁾ - **Gabriel Myers**⁽¹⁾ - **Jeffery Dahlberg**⁽²⁾ - **Joy Hollingsworth**⁽²⁾ - **Julie Sievert**⁽²⁾ - **Robert Hutmacher**⁽³⁾ - **Mary Madera**⁽⁴⁾ - **Peggy Lemaux**⁽⁴⁾ - **Kim Hixson**⁽¹⁾ - **Christer Jansson**⁽¹⁾ - **Ljiljana Paša-Tolić**⁽⁴⁾

Pacific Northwest National Laboratory, Environmental Molecular Sciences Laboratory, Richland, United States (1) - University of California, Kearney Agricultural Research & Extension Center, Parlier, United States (2) - University of California - Davis, Westside Research & Extension Center, UC Department of Plant Sciences, Five Points, United States (3) - University of California - Berkeley, Department of Plant and Microbial Biology, Berkeley, United States (4)

Keywords: top-down mass spectrometry, plant, epigenetics, post-translational modifications, molecular biology

Introduction

Sorghum is a cereal crop noted for its ability to survive water-limiting conditions. We are studying epigenetic control of drought responses in sorghum to understand how to equip similar crops to better survive climate change, and also potentially design ways to increase yield under low-water conditions. This study focuses on understanding changes in sorghum histone modifications through developmental stages and in response to drought stress in two genotypes.

Methods

Sorghum leaf tissue was collected in the field at Kearney Agricultural Research & Extension Center (Parlier, CA, USA) and frozen in liquid nitrogen on-site. Nuclei were extracted from ground leaf powder and lysed. Histones were purified from the lysate with weak cation exchange chromatography and directly analyzed with reversed phase liquid-chromatography (Waters NanoAcquity) mass spectrometry (Thermo LUMOS Orbitrap). Data were processed with TopPIC. [1]

Results

Sorghum plants from pre-flowering (PreF) and post-flowering (PostF) drought-tolerant genotypes, RTx430, and BTx642, respectively,[2] were replicated in 3 plots. All plots were watered prior to planting. The control and PostF plots were watered from the 3rd week. The PreF plots experienced 8 weeks of drought from planting to flowering. The PostF plots experienced drought after flowering in week 9 through harvest. We analyzed leaf histones from the 3rd, 9th and 10th week. LCMS experiments identified >300 histone proteoforms out of ~30 histone proteins (unique FASTA entry) per sample, including two distinct subtypes of H2A around 16 and 14 kDa, respectively. Most notably, comparative analysis revealed significant differences in the relative abundances of truncated H4 proteoforms, which showed highly elevated levels only at week 9 for leaves from control and PreF plots from BTx642, and week 10 for PreF in RTx430. The distinct responses imply that H4 clipping is related to post-flowering development and epigenetic control of drought resistance.

Conclusions

Earlier studies have suggested regulatory functions of H3/H4 clipping on nucleosome conformation and dynamics. [3] However, it is challenging to characterize histone clipping and other combinatorial

modifications using standard bottom-up proteomics methods [4–6]. Herein, we demonstrate the potential of top-down MS in profiling combinatorial histone modifications to improve our understanding of these elusive epigenetic control mechanisms.

Novel Aspect

First comprehensive top-down study of histones extracted from sorghum plants grown in the field reveals significant changes of modifications related to drought resistance.

References

1. Kou, Q., Xun, L., Liu, X. : *Bioinformatics*. 32, btw398 (2016). doi:10.1093/bioinformatics/btw398
2. Xu, L., Naylor, D., Dong, Z., Simmons, T., Pierroz, G., Hixson, K. K., Kim, Y. -M., Zink, E. M., Engbrecht, K. M., Wang, Y., Gao, C., DeGraaf, S., Madera, M. A., Sievert, J. A., Hollingsworth, J., Birdseye, D., Scheller, H. V., Hutmacher, R., Dahlberg, J., Jansson, C., Taylor, J. W., Lemaux, P. G., Coleman-Derr, D. : *Proc. Natl. Acad. Sci. U. S. A.* 115, E4284–E4293 (2018). doi:10.1073/pnas.1717308115
3. Nurse, N. P., Jimenez-Useche, I., Smith, I. T., Yuan, C. : *Biophys. J.* 104, 1081–8 (2013). doi:10.1016/j.bpj.2013.01.019
4. Tvardovskiy, A., Wrzesinski, K., Sidoli, S., Fey, S. J., Rogowska-Wrzesinska, A., Jensen, O. N. : *Mol. Cell. Proteomics*. 14, 3142–3153 (2015). doi:10.1074/mcp.M115.048975
5. Yuan, Z. -F., Arnaudo, A. M., Garcia, B. A. : *Annu. Rev. Anal. Chem.* 7, 113–128 (2014). doi:10.1146/annurev-anchem-071213-015959
6. Siuti, N., Kelleher, N. L. : *Nat. Methods*. 4, 817–821 (2007). doi:10.1038/nmeth1097

TPS-S07 - Mass spectrometric analysis of forensic science evidence

Chairs: Simona Francese, Gloriano Moneti

TP-151 / HUMAN SALIVARY GLYCOMICS FOR SALIVA-SPECIFIC SIGNATURES BY NANOPGC-LC/MS

Bum Jin Kim⁽¹⁾ - **Jin Young Park**⁽²⁾ - **Hyun Joo An**⁽²⁾

Asia-Pacific Glycomics Reference Site, Chungnam National University, Daejeon (1) - Graduate School of Analytical Science and Technology & Asia-Pacific Glycomics Reference Site, Chungnam National University, Daejeon (2)

Keywords: Saliva, Glycosylation, Mass spectrometry, *N*-glycan, Forensic

Introduction

Human saliva, which is often found at the scene of the crime, must be distinguished from other fluids in order to be recognized as legal evidence. However, conventional methods to identify saliva have lack of specificity due to cross-reactivity with other fluids. Here, we have explored saliva glycosylation using nanoLC/MS to examine saliva-specific signatures which can differentiate human saliva from other biological fluids.

Methods

Whole saliva was collected from 30 healthy donors. *N*-glycans were directly released from saliva (50µL) by PNGase F and then enriched by PGC-SPE. Saliva *N*-glycans were comprehensively analyzed by nanoLC-PGC chip/Q-TOF MS and tandem MS. To determine human saliva-specific glycans, human saliva *N*-glycans were compared with human fluids and animal saliva. Individual and environmental variations of saliva glycosylation were also assessed.

Results

Approximately 100 *N*-glycans were determined through human saliva samples. Interestingly, most glycans found in saliva were decorated with more than one fucose. In particular, highly fucosylated *N*-glycans containing at least 4 fucose residues were observed in abundance (up to 14%). Interestingly, sialylated *N*-glycans were highly abundant (>70%) in human serum. Acidic glycans with *O*-acetylation on a sialic acid were found in rat saliva as major components at a high rate while sulfation on outer HexNAc was solely observed in bovine saliva. Conclusively, highly fucosylated *N*-glycans (n=19) and significant fucosylation levels (> 50%) were determined as saliva-specific signatures for identifying human saliva from other vital fluids. These significant glycans were further investigated to validate individual variation and sample stability. Human saliva-specific glycans were consistently observed regardless of the collection date while they were gradually decreased during storage time of saliva (1, 2, 3, and 7 days) at room temperature.

Conclusions

This study demonstrated the feasibility of MS-based glycomic approach for distinguishing human saliva from other body fluids. Nineteen highly fucosylated *N*-glycans were exclusively present in human saliva and fucosylated *N*-glycans account for more than 50% quantity of total glycans found in saliva. Saliva-specific glycans and their amounts can be used as a bio-signature to differentiate human saliva from other biological fluids in forensic science.

Novel Aspect

For the first time human saliva-specific glycans were identified by mass spectrometry-based glycomic approach.

References

1. Bum Jin K., Chanyoung H., Hantae M., Joseph K., Ik-Soon J., Si-Keun L., Ki-Won P., Jong-Soon C. Hyun Joo A., *Analytical and Bioanalytical Chemistry*, 410, 45-56 (2017).
2. Sureyya O., Hyun Joo A., Ana C. V., Gun Wook P., Jae Ham K., Mark J. M., Carlito B. L., *Journal of Proteome Research*, 12, 1090-1100 (2013).
3. Hyun Joo A., Jae Han K., Bum Jin K., Jong-Soon C., Dong-Gi L., Joseph K., Verifying method for human saliva using human saliva-specific glycan, KR10-1672717 (2016), patented.

TP-152 / SINGLE-INJECTION SCREENING OF 664 FORENSIC TOXICOLOGY COMPOUNDS USING AN INNOVATIVE BENCHTOP HIGH RESOLUTION MASS SPECTROMETER

Janna Anichina⁽¹⁾ - *Oscar G. Cabrices*⁽²⁾ - *Xiang He*⁽²⁾ - *Holly McCall*⁽²⁾ - *Laura Baker*⁽²⁾ - *Alexandre Wang*⁽²⁾ - *Adrian Taylor*⁽³⁾

SCIEX, Product Applications Laboratory, Concord (1) - SCIEX, Applications, Redwood City (2) - SCIEX, Product Management, Concord (3)

Keywords: Quadrupole Time-of-flight mass spectrometry, drug screening, single injection

Introduction

Quadrupole Time-of-flight mass spectrometry (QTOF-MS) provides high-resolution, accurate-mass data for full-scan information of both precursor and product ions making it ideal for forensic screening where unknown compounds in biological samples must be identified. In this presentation, a single-injection method for screening for the enhanced identification of 664 most up-to-date forensic compounds using an innovative benchtop QTOF-MS system is described

Methods

Urine and whole blood samples were spiked with stock standard mixtures and used to determine the retention times of the 664 compounds. Samples were extracted and further diluted with mobile phase for

chromatographic separation. Mobile Phase was ammonium formate in water and formic acid in methanol, 600 $\mu\text{L}/\text{min}$ flow rate. The QTOF-MS was operated in positive electrospray mode with information dependent acquisition MS/MS methods for accurate compound identification.

Results

Separation performance was evaluated with different mobile phases (acidic and neutral), gradient conditions, and column types. Results indicate that most of isomeric compounds were fully resolved with a neutral buffer and a 10-min linear gradient. Analyte retention time (RT) was a critical element for accurate identification of all 664 compounds using this screening method. Reproducibility tests indicated that the RTs obtained from the optimized LC conditions were consistent and reproducible. RTs measured have %CVs of less than 5% for all compounds. RT inter-day reproducibility resulted in %CVs less than 5% over 3 days. RT variability in human whole blood and urine samples indicated that the %CV for 3 individual lots is less than 5%. The optimized LC conditions combined with MS/MS spectra, enabled accurate compound identification with this workflow. Retrospective analysis was also performed on the acquired data sets to screen for new compounds without having to re-inject samples, allowing to re-process data sets as new targets were discovered.

Conclusions

We have developed a QTOF-MS screening method for 664 forensic compounds, the optimized LC conditions combined with MS/MS spectra, enabled accurate compound identification. In addition, because the data was acquired in a non-targeted approach the processing method designed here for screening targeted compounds can be quickly adjusted and used for unknown compound identification using a non-targeted data processing, without having to re-inject samples.

Novel Aspect

High resolution Mass Spectrometry Analysis of 664 forensic compounds in a single injection using an innovative benchtop QTOF-MS system.

TP-153 / ULTRA-SENSITIVE FORENSIC ANALYSIS WORKFLOW OF COCAINE AND METABOLITES IN HAIR SAMPLES USING LC-MS/MS

Janna Anichina⁽¹⁾ - *Oscar G. Cabrices*⁽²⁾

SCIEX, Product Applications Laboratory, Concord (1) - SCIEX, Applications, Redwood City (2)

Keywords: triple quadrupole linear ion trap mass spectrometry, Cocaine metabolites, hair analysis

Introduction

Forensic analysis of cocaine in hair requires a sensitive and reliable analytical workflow. There are two major challenges for the detection of this compound and its metabolites in hair samples: Low concentration and matrix interferences. In this presentation, an analysis workflow combining the use of triple quadrupole linear ion trap mass spectrometry with solid phase extraction (SPE) for the picogram detection of Cocaine and Metabolites in hair is described.

Methods

Hair samples were washed, dried and cut into ~ 2 mm segments. The analytes were extracted via sonication for four hours. Supernatants were subsequently passed through SPE cartridges, dried and reconstituted in mobile phase for HILIC separation. Mobile Phase was ammonium formate in water and acetonitrile, 400 $\mu\text{L}/\text{min}$ flow rate. A SCIEX QTRAP[®] 6500+ LC-MS/MS system operating in positive electrospray mode was used for detection, with two MRM transitions for all analytes.

Results

One of the main challenges encountered throughout the study was to obtain of clean hair extracts for the reliable detection and quantitation of cocaine and metabolites from hair due to the presence of complex matrix contents. To remove these interferences, Phenomenex Strata-X SPE cartridges were successfully used to selectively isolate the drug and its metabolites. Analyte extraction recoveries were demonstrated to be greater than 80 % enabling the analytical workflow to obtain sub pg/mg Lower Limits of Quantitation (LLOQ) in hair matrix. The workflow showed excellent accuracy (>95%) and precision (< 15%), with excellent linearity resulting in R2 values of 0.9990 for all analytes. In addition to quantitation, the SCIEX QTRAP® 6500+ LC-MS/MS system enabled simultaneous identification and confirmation of illicit drugs and their metabolites through utilization of Enhanced Product Ion Scan (EPI) by acquiring full scan MS/MS data. Forensic drug identification and confirmation was achieved using automated MS/MS library searching.

Conclusions

A quantitative SPE-LC/MS/MS analysis workflow for simultaneous determination of Cocaine and metabolites at picogram levels was developed and evaluated in hair matrix. The method was demonstrated to be sensitive, precise and accurate. Utilization of the SCIEX QTRAP® 6500+ LC-MS/MS system was demonstrated to provide unique advantages in the ability to maximize selectivity when confirming and quantifying low level metabolites in hair.

Novel Aspect

An efficient analysis workflow designed with triple quadrupole linear ion trap mass spectrometry allows the picogram detection of the cocaine and metabolite in hair samples

TP-154 / FT-ICR FINGERPRINTING OF HUMAN FINGERPRINTS

Ana Cristina Assis⁽¹⁾ - *Maria Helena Florêncio*⁽²⁾ - *Carlos Cordeiro*⁽²⁾

Laboratório de Polícia Científica da Polícia Judiciária, University of Lisbon/Laboratório de FTICR e Espectrometria de Massa Estrutural, FCUL, Lisbon (1) - University of Lisbon, University of Lisbon/Laboratório de FTICR e Espectrometria de Massa Estrutural, FCUL, Lisbon (2)

Keywords: FT-ICR-MS, ESI, Lipids, Forensic Analysis, Fingerprints

Introduction

Identification of latent fingerprint composition has been carried out by several techniques. The technological development of recent years has boosted the increasingly visible use of mass spectrometry. However, so far, the chemical / biochemical fingerprint analysis has been performed by conventional mass spectrometry methods.

Methods

In the present work we characterized the lipidic profile of fingerprints on a non-porous surface by direct analysis using ESI-FT-ICR mass spectrometry and demonstrate the applicability of ESI-FT-ICR-MS with Electrospray Ionization (ESI) in the detection and characterization of lipids from human fingerprints.

Results

In a single analysis, thousands of chemical species of different origins can be identified, thus generating extremely valuable information that will lead to a chemical signature of a person at a given time. Concerning lipids, 140 lipids were identified, corresponding to 13 FA sub-classes, 6 PR sub-classes, 17 ST sub-classes, 12 PK sub-classes, 2 SP sub-classes and 1 GL sub-class, making a total of 51 sub-classes.

Conclusions

This technique uncovered a high number of lipid subclasses, both of endogenous and exogenous origin, like

what will be obtained in real cases. This characterization will allow the identification of individualizing and differentiating elements to generate forensic evidence in criminal investigations.

Novel Aspect

FT-ICR-MS analysis of human fingerprint.

TP-155 / A QUALITATIVE/QUANTITATIVE LC-QTOF-MS ASSAY FOR FORENSIC DRUG SCREENING IN URINE – FEASIBILITY STUDY AND BASIC METHOD VALIDATION

Laura Huppertz⁽¹⁾ - *Ronja Peter*⁽¹⁾ - *Michaela Schmidt*⁽¹⁾ - *Zoltan Czentnar*⁽²⁾ - *Volker Auwärter*⁽¹⁾

University of Freiburg, Institute of Forensic Medicine, Freiburg (1) - Bruker Daltonik GmbH, Application/ Applied, Bremen (2)

Keywords: HRMS, LC-QTOF-MS, qual/quant, screening, retrospective

Introduction

For forensic analysis, screening methods using LC-QTOF-MS are valuable tools due to the possibility of qualitative/quantitative and retrospective data evaluation in a single run. Here, a workflow was validated for qualitative and quantitative analysis of drugs and drugs of abuse in human urine and to assess its limitations and to prove that cut-off values are met, a validation including LOD, LOQ, linearity, accuracy, selectivity, and precision was carried out.

Methods

93 substances of forensic relevance were spiked into ten different urine samples at concentrations 1 - 50 ng/mL. Samples with internal standards were precipitated with acetonitrile. The dried residues were reconstituted and analyzed on three LC-MS systems (Bruker Elute UHPLC with impact II QTOF) in ESI pos. mode with full scan and broadband CID spectra. Data evaluation was performed with TASQ using a database with MS and chromatographic information on 2184 compounds.

Results

Identification and semi-quantification at the lowest concentration (1.0 ng/mL) was achieved for 60 % of the compounds; four compounds (paracetamol, norclomipramine, piritramid, and levetiracetam) were not detected, probably due to matrix effects and/or low ionization yields. LOD was set to the concentration at which a substance was detected in 95% of all measurements according to the following identification criteria: retention time ± 0.3 min, signal to noise ratio $> 3:1$ for all ions, $[M+nH]^{n+}$ and $[M+nH+1]^{n+}$ detected (MS), and at least two qualifier ions with minimum one being a true fragment of the molecular ion (bbCID). The precision ranged from 8 % to 30 %. Overall accuracy was good. Noteworthy, designer benzodiazepines and fentanyl derivatives were detected with extraordinarily high sensitivity. Typical 'date rape drugs' like flunitrazepam, doxylamine, and diphenhydramine showed LODs sufficient for detection of a recent uptake of these drugs.

Conclusions

All substances with legal cut-offs according to German regulations were detected below the respective concentrations (factor 5).

Extrapolating the presented urine analysis results, application to blood and hair samples seems promising and is to be evaluated/validated. Given the high frequency of new psychoactive substances emerging on web-based drug markets, this is of great interest to the forensic field due to the possibility of retrospective data evaluation.

Novel Aspect

Fully retrospective, highly sensitive qualitative/quantitative HRMS screening assay featuring semi-automated data acquisition and evaluation.

TP-156 / ANALYSIS OF SUBSTANCES RELEVANT TO § 24A (2) OF THE GERMAN ROAD TRAFFIC ACT USING LC-MSⁿ – A METHOD DEVELOPMENT AS A FIRST STEP TO AVOID IMMUNOASSAY PRE-SCREENING.

Jürgen Kempf⁽¹⁾ - *Volker Auwärter*⁽¹⁾ - *Franziska Ehrhardt*⁽²⁾ - *Zoltan Czertnar*⁽³⁾

University of Freiburg, Institute of Forensic Medicine, Freiburg (1) - Offenburg University, Applied Sciences, Offenburg (2) - Bruker Daltonik GmbH, Application/Applied, Bremen (3)

Keywords: § 24a (2) road traffic act, method development, LC-ion trap-MS

Introduction

Cases of suspected driving under the influence of drugs (DUID) according to § 24a (2) of the German road traffic act (RTA) usually pre-screened by immunoassay (IA) and positive results are confirmed by LC-MS/MS since neither the qualitative nor the quantitative information from IAs is admissible in court. We developed a fast and automated LC-MSⁿ method, combining the ease-of-use of IAs with unambiguous compound identification of an MS analysis.

Methods

Step one of this project covers a general feasibility study including two different data acquisition modes (AutoMSⁿ and sMRM) for qualitative and/or quantitative analysis. It contains a sample preparation of 0.5 ml serum using LLE with chloroform/isopropanol, applicable for all relevant analytes except THC. Analysis is carried out on an Acclaim C18 column using a UltiMate3000 LC-system coupled to an amaZon speed ion trap within 5 minutes runtime.

Results

During this study two data acquisition methods were evaluated and implemented: The Auto-MSⁿ-mode facilitates fast and easy qualitative analysis with automated result reporting whereas the sMRM-mode offers a first quantitative estimate of the analyte concentration, simplifying subsequent confirmation analysis. Both methods include the § 24a (2) RTA relevant compounds amphetamine, methamphetamine, ecstasy, heroin, morphine, cocaine and benzoylecgonine as well as codeine and 6-acetylmorphine and have been validated according to the GTFCh guidelines. Selectivity and LODs were verified for both methods. For the sMRM-method linearity, matrix effects, stability and accuracy were also determined. Measurement of 60 real case samples showed that neither the AutoMSⁿ- nor the sMRM-method generated false-positive or false-negative results.

Conclusions

Although sample preparation is still carried out manually at this point, the developed LC-MSⁿ approach would be a suitable replacement for IA testing in DUID cases according to § 24a (2) RTA. As a next step of this project, an online SPE is currently evaluated and will be implemented to enhance the degree of automation of this approach.

Novel Aspect

Detection of § 24a (2) RTA relevant compounds using a fast and robust LC-MSⁿ method to avoid immunoassay pre-screening in routine DUID analysis.

TP-158 / APPLICATION OF TOF-SIMS TECHNIQUE TO DETECTION AND IMAGING OF FINGERPRINTS CONTAMINATED BY CANNABIS

Elzbieta Mackiewicz⁽¹⁾ - *Jacek Rogowski*⁽¹⁾ - *Andrzej Parczewski*⁽²⁾ - *Małgorzata Iwona Szyrkowska*⁽¹⁾

Lodz University of Technology, Faculty of Chemistry, Lodz, Poland (1) - Jagiellonian University, Faculty of Chemistry, Krakow, Poland (2)

Keywords: fingerprints; imaging; TOF-SIMS; drug analysis; development of latent fingerprints

Introduction

Fingerprints are the most important form of evidence in the human identification [1]. Moreover, they may be composed of a number of different contaminants, for example psychoactive substances[2,3]. During the last few years, time-of-flight secondary ion mass spectrometry (TOF-SIMS) has become a very sensitive analytical technique in surface science [2,3]. This technique can be used in chemistry and forensic science for detection of various molecules (drugs).

Methods

The mass spectra and images of the studied samples were recorded using secondary ion mass spectrometer TOF-SIMS IV (ION-TOF GmbH, Münster, Germany), which is equipped with a high mass resolution time-of-flight analyser of a reflectron type. During analysis, Bi³⁺ primary ion gun was used, working at 25 keV and an average ion current was 0.5 pA. The analysis areas were 100 × 100 μm² and 500 × 500 μm² for obtaining secondary ion mass spectra and images, respectively.

Results

Time-of-flight secondary ion mass spectrometry (TOF-SIMS) was applied in the forensic research for detecting and imaging of cannabis in the fingerprints. The studied fingerprints were: 1. imprinted directly on the lift tape by a contaminated finger with cannabis; 2. left on a glass surface and developed with two types of latent fingerprints powders; 3. transferred using the lift tape.

In the positive and negative ion spectra of the studied samples many peaks derived from cannabis can be detected: cannabidiol, cannabidiol, tetrahydrocannabinol, cannabivarin, cannabidiolic acid. The ions characteristic of the powder, sweat and marijuana are also visible in the spectra of fingerprints revealed with latent fingerprints powder and transferred with lift tapes.

Conclusions

The TOF-SIMS technique appears to be a very effective and useful tool in forensic investigations, which enables the identification of drug residues left on fingerprints imprinted on selected bases. Moreover, it has been shown that revealing fingerprints with the black fingerprint powder is not disadvantageous for the TOF-SIMS application in the fingerprints analysis.

Novel Aspect

There is a possibility of detecting foreign substances on fingerprints, such as drugs, after revealing fingermarks with the latent fingerprint powder.

References

1. Hamilton I. Encyclopedia of Forensic Sciences 346 (2013).
2. Szyrkowska M. I., Czerski K., Rogowski J., Paryjczak T., Parczewski A., Surface and Interface Analysis 42, 393 (2010).
3. Szyrkowska M. I., Czerski K., Rogowski J., Paryjczak T., Parczewski A., Forensic Science International 184 e24 (2009).

TP-159 / VALIDATION OF BLOOD PROTEOMIC SIGNATURES FOR THE FORENSIC DETERMINATION OF BLOOD AND ITS PROVENANCE

Cameron Heaton⁽¹⁾ - ***Lisa Deininger***⁽¹⁾ - ***Ekta Patel***⁽¹⁾ - ***Malcolm. R. Clench***⁽¹⁾ - ***Vaughn Sears***⁽²⁾ - ***Chris Sammon***⁽³⁾ - ***Laura Cole***⁽¹⁾ - ***Simona Francese***⁽¹⁾

Sheffield Hallam University, Biomolecular Research Centre, Sheffield (1) - Home Office, Centre for Applied Science and Technology, St Albans (2) - Sheffield Hallam University, Materials and Engineering Research Institute, Sheffield (3)

Keywords: Blood, Forensics, Proteomics, MALDI

Introduction

Analysis of blood at a crime scene can inform an investigation but current enhancement techniques can provide false positives. Bottom up proteomic MS can detect blood and indicate its origin [1]. Blood specific proteins can be visualised directly on the ridges of a fingerprint [2]. In this study, human and animal blood, biofluids, and substances that elicit a positive response to presumptive tests have been examined pre- and post-enhancement. Blood identification was based on proteotypic peptides. Manual and automated PMF methods were trialled to identify the samples.

Methods

Blood stains and fingerprints on glass slides were received from Elite Forensics[®]. Stains were digested with 250 µg/ml trypsin after swabbing and extraction, followed by spot analysis via MALDI MS Profiling. Fingerprints were digested in situ using a SunCollect auto-spraying system, incubated for 1 hour and then sprayed with matrix (αCHCA) using the same system. MALDI analyses were carried out on a Synapt G2 HDMS. Images were processed using Waters HDI software.

Results

Preliminary results show a 100% correct identification rate when human blood samples were analysed blindly in enhanced and non-enhanced form. This classification rate is based on 1/3 of the samples analysed so far. The prior use of blood enhancement techniques has not interfered with the retrieval of proteotypic peptides. Other biofluids and substances such as egg yolk were correctly classified as non-blood. Manual PMF is based on specific criteria including ion signal exclusion from the mass lists. Whilst manual PMF is highly successful in sample classification, the use of the automatic PMF software correctly attributes human blood with a low probability score. Animal blood presents additional difficulties. Prior investigations have correctly classified animal blood. In this study, some animal blood samples were classified as non-blood, likely to be as a result of blood sampled from meat packets rather than venous blood, which has led to the acquisition of spectra with low signal to noise impacting on the PMF, both in manual and automatic mode.

Conclusions

This first stage validation study has developed methods to correctly classify human blood and non-blood samples. Whilst correct classification of animal blood is possible, determination of provenance has been compromised. Further research into improvement of MS acquisition methods as well as of the treatment of the data is in progress. Venous animal blood (instead of from meat packets) will also be trialled and results compared to assess robustness of the methods.

Novel Aspect

Proteomics MALDI MS methods were validated for use in police casework. Robustness of PMF software using a peptide library was assessed against manual PMF of proteotypic peptides.

References

1. Patel E., Cicatiello P., Deininger L., Clench M. R., Marino G., Giardina P., Langenburg G., West A., Marshall P., Sears V., Francese S., *Analyst*, 141, 191 (2016)
2. Deininger L., Patel E., Clench M. R., Sears V., Sammon C., Francese S., *Proteomics*, 16, 1707 (2016)

TP-160 / DEVELOPMENT AND VALIDATION OF A GAS-CHROMATOGRAPHY MASS SPECTROMETRY METHOD FOR THE ANALYSIS OF TESTOSTERONE ESTERS IN BIOLOGICAL MATRICES

Michele Iannone⁽¹⁾ - ***Claudia Lattanzi***⁽¹⁾ - ***Francesco Botrè***⁽²⁾ - ***Xavier de la Torre***⁽¹⁾

Laboratorio Antidoping FMSI Rome, Dipartimento di Chimica e Tecnologia del Farmaco, Sapienza University

of Rome, Rome (1) - Laboratorio Antidoping FMSI Rome, Dipartimento di Medicina Sperimentale, Sapienza University of Rome, Rome (2)

Keywords: Testosterone esters, doping analysis, gas-chromatography, mass spectrometry.

Introduction

The detection of pseudoendogenous steroid doping is currently based on the evaluation of the urinary steroid profile as described in WADA Technical Document TD2016EAAS [1]. Despite this some physiological and/or pathological conditions it could compromise the application of these criteria [2-3]. For this reason the direct detection of intact synthetic esters of testosterone in plasma/serum gives unequivocal proof of exogenous testosterone administration.

Methods

After liquid-liquid extraction at pH 9 (carbonate buffer 20% w/V) with n-hexane/ethylacetate (70/30 V/V) [4] all plasma/serum samples were analyzed by gas-chromatography coupled to tandem mass spectrometry (GC-MS/MS) with electronic impact (EI) ionization and multiple reaction monitoring (MRM) acquisition mode. Testosterone esters were revealed as trimethylsilyl derivatives.

Results

In the first step of our work a qualitative validation was carried out for all testosterone esters studied (acetate, propionate, valerate, isocaproate, caproate, benzoate, hexahydrobenzoate, enanthate, cypionate, octanoate, decanoate, phenylpropionate, undecanoate, laurate). Considered parameters are specificity/selectivity and limit of detection. According to the WADA rules [5] were also investigated the mass spectrometric and the chromatographic criteria for the identification of analytes for doping control purposes. The criteria were satisfied for twelve compounds and the method is specific and selective for all of these.

In a second step a quantitative validation was carried out for five compounds present in the most marketed pharmaceutical preparations probably abused by athletes. Considered parameters are specificity/selectivity, limit of detection and quantification, recovery, linearity, repeatability and accuracy. The selected esters are identified and quantified in plasma/serum with an error less than 15 %.

Conclusions

The determination of testosterone esters in plasma/serum by GC-MS/MS should allow a confirmation of the administration of exogenous testosterone for doping purpose. The presented method was qualitative validated for twelve esters and quantitative validated for five esters and can be used to detect the presence in blood of those testosterone esters currently present in pharmaceutical preparations that can be abused by athletes.

Novel Aspect

Number of esters detected in each sample and possibility to confirm the composition of seized pharmaceutical preparation by responsible authorities.

References

1. WADA Technical Document TD2016EAAS. Endogenous Anabolic Androgenic Steroids. Measurement and Reporting. Montreal, September 2015.
2. Mareck U., Geyer H., Opfermann G., Thevis M., Schanzer W. *J Mass Spectrom*, 43, 877-891 (2008).
3. Kuuranne T., Saugy M., Baume N. *Br J Sports Med*, 48, 848-855 (2014).
4. De la Torre X., Segura J., Poletti A., Montagna M. *J Mass Spectrom*, 30, 1393-1404 (1995).
5. WADA Technical Document TD2015IDCR. Minimum Criteria for Chromatographic-Mass Spectrometric Confirmation of the Identity of Analytes for Doping Control Purposes. Montreal, September 2014.

TP-161 / ELEMENTAL ANALYSIS OF INK BY LASER ABLATION INDUCTIVELY COUPLED PLASMA MASS SPECTROMETRY (LA-ICP-MS)

Aleksandra Pawlaczyk - Karolina Załęcka - Małgorzata Iwona Szyrkowska

Lodz University of Technology, Faculty of Chemistry, Institute of General and Ecological Chemistry, Lodz

Keywords: forensic analysis, chemometric methods, inks, questioned documents, LA-ICP-MS

Introduction

In real cases the analysis of questioned documents is undertaken with mostly non-invasive methods of examination of paper base & covering materials. The methods applied should provide reliable and sufficient information to address the following questions like whether the document is genuine, if it has been altered (e.g. erased, bleached, faded) or reproduced, if the covering materials (toners, inks) appearing on each page of document have the common source [1].

Methods

ICP Mass Spectrometer with Time-of-Flight Analyzer coupled with Laser Ablation system was applied for inorganic chemical elements characterization of different inks. In this study 13 random, promotional pens were collected and used to leave a straight few cm lines on a white piece of paper. The LA-ICP-MS parameters were optimized in order to achieve satisfactory intensity of the ink signal for selected m/z , leaving the paper base almost untouched.

Results

Inks are complex mixtures of different colorants, vehicles, and additives, which depending on the company, color and type are adjusted in composition to produce the desired writing characteristics. In order to avoid loss or destruction of a potential forensic evidence hardly destructive methods like LA-ICP-MS are desirable. This technique with high discrimination and classification power can be successfully applied to differentiate inks on the same document as well as to compare inks on other documents. In this work the mass spectra profiles for each randomly selected & commercially available on the market inks (blue and black) were recorded and compared. They all presented specific elemental features which made them distinctly different from each other. The gathered information about elemental composition differences was supported by the feedback from chemometric methods. In the next step of analysis the "unknown" sets of 3 chosen samples using studied pens were prepared, tested, discriminated and classified according to created model.

Conclusions

The gathered data originating from collected mass spectra coupled with multivariate statistical methods present a useful tool for discrimination of studied unknown inks with insubstantial sample destruction and without any methodological calibrations, which all being crucial from a forensic practice point of view. The same methodology can be applied to analyze real handwritten questioned documents since there is a justified need to use non-destructive methods.

Novel Aspect

There is a constant need to create method with high discrimination power which can deliver objective results offering the differentiation among brands & within the same brand.

References

1. Chęsy P., Pawlaczyk A., Szyrkowska M. I., Parczewski A., *Problems of Forensic Sciences*, 108, 605 (2016).
2. Orellane F. A., Gálvez C. G., Roldán M. T., García-Ruiz C., *Trends in Analytical Chemistry*, 42, 1 (2013).

TP-162 / DETECTION OF RICIN AND ABRIN TOXINS IN COMPLEX MATRICES BY TWO NOVEL GALACTOSE ENRICHMENT METHODS AND HIGH RESOLUTION LC-MS/MS

Anne Puustinen - Marja-Leena Rapinoja - Maaret Karjalainen - Paula Vanninen

Verifin, Dept of Chemistry, University of Helsinki, Helsinki, Finland

Keywords: ricin, abrin, galactose affinity, LC-MS/MS

Introduction

Ricin and abrin are highly poisonous plant toxins belonging to the type 2 ribosome-inactivating protein (RIP-II) family. These proteins are composed of A- and B-chains linked by a single disulphide bond. The A-chain confers cellular toxicity, while cellular entry is through the galactose binding properties of the B-chain, which can be utilized for isolation of these toxins. There is a demand for improving toxin detection capabilities due to the terrorism threat posed by them.

Methods

Immobilized D-galactose gel (Thermo Scientific) in self-packed microspin columns and home-made galactose-coated 96-well microplates were utilized for enrichment of ricin and abrin from different kind of samples. Peptides from the trypsin digested protein samples were identified by high resolution mass spectrometry analysis (EASY nanoLC linked to Fusion Orbitrap, Thermo Scientific). Released adenine in the depurination activity assay was measured with LC-MS/MS (UPLC linked to Xevo TQD, Waters).

Results

Two rapid galactose affinity enrichment methods were developed for ricin and abrin analysis from complex matrices. Miniaturized immobilized D-galactose gel affinity columns were shown to be suitable for trace enrichment and high sample throughput. This method can also be adapted for preparative purification of RIP-II toxins. Also the solid-phase enrichment by galactose array worked for trace enrichment and can be automated with microplate automation instruments.

Ricin and abrin are fairly resistant to proteolysis. Protein unfolding was done with an acid-labile surfactant - urea mixture before digestion with trypsin. Identities of the resulting peptides were unambiguously confirmed by accurate mass LC-MS/MS product ion scans of peptides unique to each of the toxins. Quantification was accomplished with isotope labeled AQUA peptides for both proteins and the achieved sensitivity and detection limits were comparable to those of immunocapture enrichment methods.

Conclusions

Both galactosyl ligand methods were shown to have the capacity to enrich trace amounts of ricin and abrin from various matrix extracts. Together with functional assay results, the high resolution LC-MS/MS data provides unambiguous identification of ricin and abrin in the samples. The solid-phase enrichment by galactose array can be adapted for high-throughput screening.

Novel Aspect

High throughput galactose affinity methods for ricin and abrin with detection limits similar to the ones reported for other enrichment methods were developed.

TP-163 / LC-HRMS VS NMR FOR THE IDENTIFICATION OF KETOCAINE IN SEIZED MATERIAL

***Anna Maria Caputo*⁽¹⁾ - *Serena Detti*⁽¹⁾ - *Sabino Napoletano*⁽¹⁾ - *Giovanni De Sangro*⁽¹⁾ - *Rossana Scarpone*⁽²⁾ - *Fabio Sciubba*⁽³⁾ - *Flaminia Vincenti*⁽³⁾ - *Manuel Sergi*⁽⁴⁾ - *Dario Compagnone*⁽⁴⁾**

Ministry of Interior, Dep. of P. S., Forensic Police Service of Italian National Police, Drugs of Abuse Section,

Via Tuscolana 1548, Rome (1) - Istituto Zooprofilattico Sperimentale dell'Abruzzo e del Molise 'G. Caporale', Bromatology, Teramo (2) - Sapienza University, Chemistry, Rome (3) - University of Teramo, Faculty of Bioscience, Teramo (4)

Keywords: Forensic analysis; HRMS; NMR; Ketocaine

Introduction

The timeliness of analytical assessments is fundamental for Police forces, initially relying on rapid drug-screening tests. Scientific Police Service of Italian National Police adopts Young's reagent for the detection of cocaine; however, there are numerous examples of false positives that require attention to the interpretation of test result and, in some cases, confirmatory analysis must be performed. In this study an interesting case is presented.

Methods

The analytical protocol adopted for the identification of the suspect compound involved both NMR and HMRS. The sample was divided into two aliquots: one was analysed by means of LC-HRMS with Orbitrap technology and the other one was subjected to analysis by NMR. The two techniques have not been used in a complementary but competitive way: namely the analysis have been performed independently.

This approach involved both the comparison of the molecular ion and the fragmentation pattern with the database spectra by differential analysis, providing stronger identification.

Standard ^1H and ^{13}C NMR spectra were coupled with two-dimensional experiments: ^1H - ^1H TOCSY, ^1H - ^1H COSY, ^1H - ^{13}C HSQC and ^1H - ^{13}C HMBC.

Results

LC-HRMS identified Ketocaine by involved both the comparison of the molecular ion and the fragmentation pattern by differential analysis, from spectra database Compound Discoverer 2.0, using ChemSpider as source. Forensic identification require something more (more identification points), so the study of the molecule and its fragmentation was also performed. The match of theoretical fragmentation pattern with the experimental one confirmed the identification.

NMR experiments carried out in different acquisition modes (^1H - ^1H TOCSY, ^1H - ^1H COSY, ^1H - ^{13}C HSQC and ^1H - ^{13}C HMBC.) led to a suitable amount of data for a correct identification of the unknown compound: the interpretation of spectra and the signals assignment led to identification of ketocaine.

Conclusions

In this case, both NMR and HRMS proved themselves as self-sufficient techniques in the resolution of the analytical-forensic problem: in particular, HRMS, with the advantage of the speed of analysis and interpretation of data thanks to dedicated software, proves of undoubted utility for forensic investigations.

Novel Aspect

The suitability of HRMS as a self-sufficient technique in the resolution of the analytic-forensic problem was demonstrated for the identification of unknown compound.

TP-164 / TARGETED AND UNTARGETED ANALYSES BY MEANS OF LC-MS TO CHARACTERIZE THE ILLICIT SUBSTANCES ON SEIZED MATERIAL OF VEGETAL ORIGIN

Anna Maria Caputo⁽¹⁾ - *Serena Detti*⁽¹⁾ - *Sabino Napoletano*⁽¹⁾ - *Valter Fraioli*⁽¹⁾ - *Alessandro Possi*⁽²⁾ - *Alessandro Zompanti*⁽²⁾ - *Rossana Scarpone*⁽³⁾ - *Francesca Di Ottavio*⁽⁴⁾ - *Manuel Sergi*⁽⁴⁾ - *Flaminia Vincenti*⁽⁵⁾ - *Louis Felix Nothias*⁽⁶⁾ - *Pieter Dorrestein*⁽⁶⁾

Ministry of Interior, Dep. of P. S., Forensic Police Service of Italian National Police, Drugs of Abuse Section, Rome (1) - Regional Forensic Office, Italian National Police, Florence (2) - Istituto Zooprofilattico

Sperimentale dell'Abruzzo e del Molise 'G. Caporale', Bromatology, Teramo (3) - University of Teramo, Faculty of Bioscience, Teramo (4) - Sapienza University of Rome, Dept. of Chemistry, Rome (5) - University of California, Collaborative Mass Spectrometry Innovation Center, San Diego (6)

Keywords: illicit drugs, seized material, LC-MS, Molecular Network, in-silico annotation

Introduction

In the field of investigations in drugs of abuse, one of the objectives of the police forces is the characterization of the spread and the consumption of psychotropic substances: this help to trigger adequate monitoring and prevention measures[1]. Through different techniques, it is possible to identify the components present in seized sample in order to characterize the illicit substances.

Methods

A combination of targeted and untargeted analyses were performed on seized material of vegetal origin. Screening test and GC analyses showed an unusual combination of compounds. UHPLC-MS/MS analyses were performed in IDA and in MRM. LC-HRMS analyses followed by elaboration of raw data files by means of GNPS web platform and in-silico annotations in order to predict the structure of all illicit compounds.

Results

On the basis of the preliminary analyses there was the suspect of different tryptamine compounds in the seized samples so a semi-untargeted analysis by UHPLC-MS/MS were performed in IDA, showing different tryptamines that were confirmed in MRM(*i.e.* DMT, NMT, 2-MTHBC), together with some cannabinoids, including THC, in one sample. This issue stimulated further investigations.

With the aim to obtain a complete overview of all features, untargeted HPLC-HRMS analysis was performed. Molecular Network was generated to perform cluster of similar MS/MS spectra and to propagate spectral annotation of seized drug samples. Data shows that different compounds, belonging to tryptamines, were confirmed by spectral library search annotation on GNPS. Other 8 tryptamines were assigned comparing the results obtained with MRM analysis and similar pattern of fragmentation using in-silico prediction annotation.

Conclusions

Through different analytical techniques, it was possible to identify the main vegetable components present in samples. Combination of targeted/untargeted analysis and molecular network proves to be a useful tool for forensic investigations. Untargeted analysis and in-silico annotation allowed detecting 8 compounds undetected by MRM analysis and confirming the analytes detected in target analysis.

Novel Aspect

Combination of targeted/untargeted analysis, with suitable data managing, resulted necessary to characterize seized sample for subsequent evaluations.

References

1. European Monitoring Centre for Drugs and Drug Addiction (2017), Health and social responses to drug problems: a European guide, Publications Office of the European Union, Luxembourg.

TP-165 / MODAFIENDZ METABOLITES IDENTIFICATION BY MEANS OF LC-HRMS: FROM IN SILICO PREDICTION TO IN VITRO CONFIRMATION WITH RAT HEPATOCYTES

Flaminia Vincenti⁽¹⁾ - **Federico Fanti**⁽²⁾ - **Adolfo Gregori**⁽³⁾ - **Anna Rita Togna**⁽⁴⁾ - **Dario Compagnone**⁽²⁾ - **Manuel Sergi**⁽²⁾

Sapienza - University of Rome, Department of Chemistry, Rome (1) - University of Teramo, Faculty of Bioscience and Technology for Food, Agriculture and Environment, Teramo (2) - Arma dei Carabinieri,

Department of Scientific Investigation (RIS), Rome (3) - Sapienza University of Rome, Department of Physiology and Pharmacology Vittorio Erspamer, Roma (4)

Keywords: NPS, Nootropics substances, Metabolites, LC-HRMS

Introduction

New psychoactive substances (NPS) represent a tough challenge, since analytical standards are not available or included in libraries [1]. When metabolic pathway is unknown, determine target analytes is critical, so metabolites identification become fundamental for clinical and forensic investigations [2]. The aim of this work was to identify the metabolites of Modafiendz, a synthetic nootropic [3].

Methods

In-silico experiments were carried out by means of MetaPrint2DReact program, used to predict possible metabolites and cultures of rat hepatocytes were used for in vitro experiments. In order to obtain information about chemical structure of metabolic pathway, LC-HRMS experiments were carried out by means of Orbitrap analyzer, while analytes separation was carried out using a PFP-C18 Column.

Results

Once the smile code was introduced in the in silico program, a total of 15 metabolites were predicted: 7 of phase I and 8 of phase II. The incubation with hepatocytes allowed to mimic the physiological environment of the liver [4], in which the substances would be found, obtaining a realistic picture of the metabolic pathway for Modafiendz. The results obtained in silico allowed to hypothesize the metabolites structures and masses to be searched in HRMS. As a positive control of the real metabolic activity of the cells, Modafinil and its metabolites were used: Modafinil Acid and Modafinil Sulfone have been identified. In the samples obtained from in vitro experiments 6 metabolites were identified, 2 belonging to phase I and 4 belonging to phase II. In order to characterize the observed metabolites, analyzes were performed in PRM at different collision energies (NCE).

Conclusions

In this work the chemical structures of the Modafiendz metabolites were identified using rat hepatocyte followed by UHPLC-HRMS analysis. The most produced metabolites are dehalogenated derivatives obtained from an initial oxidation reaction, as regards the phase I metabolites and the metabolites conjugated with glucuronic acid or methylate, deriving from the phase II reactions.

Novel Aspect

For the first time metabolic pathway of Modafiendz was identified. This information is essential to develop methods to confirm the use of this NPS [5].

References

1. Gerostamoulos D., Elliot S., Walls H. C., Peters F. T., Lynch M., et al. *Journal of Analytical Toxicology*, 40, 318-320 (2016).
2. Mardal M., Kinyua J., Ramin P., Miserez B., et al. *Drug Testing and Analysis*, 9, 106-114 (2017).
3. Dowling G., Kavanagh P. V., Talbot B., O'Brien J., et al. *Drug Testing and Analysis*, 9, 518-528 (2017).
4. Montesano C., Vannutelli G., Fanti f., Vincenti F., et al. *Journal of Analytical Toxicology*, 41, 688-697 (2017).
5. Conceiro M., Castaneto M., Kronstand R., Huestis M. A., *Journal of Chromatography A*, 1397, 32-42 (2015).

TP-166 / IDENTIFICATION OF FINE PLASTIC MATERIALS BY THERMAL DESORPTION AND PYROLYSIS COMBINED WITH DART-MS (TDP/DART-MS)

Chikako Takei - Kenichi Yoshizawa

BioChromato, Inc., Application support, Fujisawa, Japan

Keywords: Criminalistics, foreign substances analysis, TDP/DART-MS

Introduction

Evidentiary materials in criminal come in many forms such as fine segments of fiber or bits of resin. Adhesive collection sheets are often used, but even if a fine particle can be collected, the instrumental analysis is affected by the adhesive compound adhered to a fine particle slightly, so the difficulty of distinguishing task can also further increase. The purpose of this work is to identify fine plastic materials without confounding influence from the adhesive tape.

Methods

The samples of this study were “resins sticking to adhesive tape”. Thermal desorption and pyrolysis(TDP) combined with DART-MS was composed of Q-TOF mass spec equipped with DART ion source and ionRocket TDP device. The samples were cut to ~1 mm² and put into the sample pot. TDP temperature condition was room temperature to 600°C, a rate of 100°C/min. To structurally and compositionally characterize samples, Kendrick Mass Defect (KMD) analysis was used.

Results

Nylon-6, polylactic acid(PLA) and polyethylene glycol(PEG) were used as fine plastic materials and scotch tape was used as the adhesive tape. In the mass spectra at 300–400°C (the pyrolysis reaction starting temperature), the complicated spectra of the samples show that the fine plastics were mixed with scotch tape, making it difficult to distinguish between them. From the KMD plots, the difference between the scotch tape and fine plastics were clearly identified. For the nylon-6 fine plastic, a plot of 113 Da intervals derived from its repeating structure was plotted at a different position from that of the peaks originating from the scotch tape. In a similar manner, plots with a 72 Da interval for PLA and a 44 Da interval for PEG were observed at different positions from that of the scotch tape. By setting the Kendrick mass unit based on the mass of the repeating structure, compounds with common repeat units could be aligned in the horizontal direction on the KMD plot, whereas compounds with different structures were shifted.

Conclusions

TDP/DART-MS enables the direct and rapid analysis of plastic samples. A combination of TDP/DART-MS and KMD analysis enables the identification of fine plastics without the need to detach foreign matter collected from adhesive tape. Thus, the combination of TDP/DART-MS and KMD analysis can be used for material identification of evidence in criminal cases.

Novel Aspect

A combination of TDP/DART-MS and KMD analysis is useful for material identification of evidence in criminal cases.

TP-167 / TOF-SIMS DETERMINATION OF DETECTION LIMIT OF AMPHETAMINE IN ‘STREET DRUGS’ AND POLLUTED FINGERPRINTS

Magdalena Gajek⁽¹⁾ - *Elzbieta Mackiewicz*⁽¹⁾ - *Jacek Rogowski*⁽¹⁾ - *Andrzej Parczewski*⁽²⁾ - *Małgorzata Iwona Szyrkowska*⁽¹⁾

Lodz University of Technology, Faculty of Chemistry, Lodz, Poland (1) - Jagiellonian University, Faculty of Chemistry, Krakow, Poland (2)

Keywords: Secondary ion mass spectrometry, fingerprints, drugs, forensic research, image technique

Introduction

ToF-SIMS technique is a very effective and powerful method to examine the surface composition of different materials, including those of key importance in the forensic chemistry i.a. drugs [1-2]. Currently, drugs sold in the streets contain a variety of different 'fillers'. It is important to develop a method which has discrimination power and will allow simultaneous identification of all substances present in the mixture, even in very small quantities.

Methods

The aim of this study is focused on the evaluation of the TOF-SIMS technique as a potential tool applied in the forensic analysis of fingerprints. Studied fingerprints were contaminated with trace amounts of amphetamine and other substances which are typically used as 'fillers'. These compounds do not occur naturally on fingertips and can present a significant value as evidence linking the suspect with the crime scene.

Results

In this work, clean fingerprints as well as contaminated by drug were investigated. Moreover, in order to determine the limit of detection of amphetamine, the analysis of drug tablets containing various amounts of mixtures of amphetamine and caffeine (filler) in different proportions was performed. This assessed LOD was 0,625%. Also fragmentative ions derived from amphetamine and caffeine were examined, which undoubtedly indicted their potential origin. Additionally, the ToF-SIMS technique allowed the visualization of the fingerprint contaminated by a trace amount of powdered drug.

Conclusions

The results of this study has demonstrated that the ToF-SIMS technique is a powerful tool which enables the identification of drug traces left on a contaminated fingerprint. This method allows the determination of very small amounts of the drug. Moreover, ToF-SIMS technique enables not only the identification of the substance forming the fingermark but also the visualization of very small areas of the fingerprint, including the course of the fingerprints.

Novel Aspect

Novel aspects are: designing a method enabling the assessment of LOD of amphetamine in home made tablets and its identification and visualization in the fingerprints course.

References

1. Szykowska M. I., Czerski K., Rogowski J., Paryjczak T., Parczewski J., Forensic Sci. Int., 184, 24-26, 2009
2. Szykowska M. I., Czerski K., Rogowski J., Paryjczak T., Parczewski J., Surface and Interface Analysis, 42, 393-397, 2010

TPS-S08 - Ion activation and dissociations

Chairs: Vicki Wysocki, Gianluca Giorgi

TP-168 / TANDEM MASS SPECTROMETRY TO ASSIST THE DESIGN OF INFORMATION-CONTAINING SYNTHETIC POLYMERS: AN EFFICIENT APPROACH TO FACILITATE THEIR SEQUENCING AND INCREASE THEIR STORAGE DENSITY

Jean-Arthur Amalian⁽¹⁾ - ***Gianni Cavallo***⁽²⁾ - ***Salomé Poyer***⁽¹⁾ - ***Jean-François Lutz***⁽²⁾ - ***Laurence Charles***⁽¹⁾

Institut de Chimie Radicalaire, Aix-Marseille Université, Marseille (1) - Chimie Moléculaire de Précision,

Keywords: Tandem mass spectrometry; sequencing; sequence-controlled polymers; molecular design; information-containing macromolecules

Introduction

MS/MS sequencing was found to be the most efficient methodology to retrieve digital information stored in sequence-defined synthetic polymers [1]. In these macromolecules, two comonomers of different mass are defined as the 0- and 1-bit of the ASCII alphabet and placed at defined location throughout the chain to build binary messages [2]. Since the comonomer structure dictates the dissociation behavior, it can be optimized to achieve best MS/MS readability.

Methods

Poly(alkoxyamine phosphodiester)s were synthesized using an orthogonal strategy that employs successively phosphoramidite and radical-radical coupling steps [3]. Methanolic solutions of polymers were supplemented with ammonium acetate prior to electrospray ionization in the negative ion mode. High resolution MS and MS² experiments were performed on a Synapt G2 Q-TOF mass spectrometer (Waters, Manchester, UK).

Results

Sequence-defined poly(phosphodiester)s can store large amounts of information as they can be produced in an automated manner as long monodisperse chains [4]. However, extensive fragmentation of all phosphate bonds in each monomer leads to very busy MS/MS spectra that do not systematically ensure full sequence coverage. To simplify CID data and hence facilitate this sequencing step, a weak alkoxyamine bond was introduced between all monomers: as expected in such poly(alkoxyamine phosphodiester)s, fragmentations specifically occurred at C–ON bonds, making all phosphate groups MS/MS silent [5]. Additional structural optimizations included tailoring the distance between phosphate moieties in order to promote their simultaneous deprotonation, which allowed production of fragments which charge state increases with their polymerization degree. Finally, storage density of these chains was further increased by introducing two instead of one coding moieties in each monomer without modifying their optimized dissociation behavior [6].

Conclusions

In contrast to biomolecules which structure is set by biology, sequence-defined synthetic polymers can be specifically designed to control their CID behavior and ensure full sequence coverage. So simplified dissociation rules can then be easily implemented in the recently-introduced MS-DECODER software [7], allowing efficient MS/MS data analysis and message decoding in the millisecond scale.

Novel Aspect

MS/MS assisted design of sequence-controlled synthetic polymers to improve readability of molecularly encoded information and increase their storage density.

References

1. Charles L., Laure C., Lutz J. -F., Roy R. K., *Macromolecules*, 48, 4319-4328 (2015)
2. Roy R. K., Meszynka A., Laure C., Charles L., Verchin C., Lutz J. -F., *Nat. Commun.*, 6, 7237 (2015)
3. Cavallo G., Al Ouahabi A., Oswald L., Charles L., Lutz J. F., *J. Am. Chem. Soc.*, 138, 9417-9420 (2016)
4. Al Ouahabi A., Kotera M., Charles L., Lutz J. -F., *ACS Macro Lett.*, 4, 1077-1080 (2015)
5. Amalian J. -A., Al Ouahabi A., Cavallo G., König N. F., Poyer S., Lutz J. -F., Charles L., *J. Mass Spectrom.*, 52, 788-798 (2017)
6. Cavallo G., Poyer S., Amalian J. -A., Dufour F., Burel A., Carapito C., Charles L., Lutz J. -F., *Angew. Chem.*, accepted.
7. Burel A., Carapito C., Lutz J. F., Charles L., *Macromolecules*, 50, 8290-8296 (2017)

TP-169 / FACILE ELECTRON CAPTURE DISSOCIATION ON Q-TOFS AND ORBITRAPs FOR TOP- AND MIDDLE-DOWN PROTEOMICS AT CHROMATOGRAPHIC TIME SCALES

Joseph Beckman - Valery Voinov

e-MSion, Inc, Oregon State University, Corvallis

Keywords: ECD, native, de novo, sequencing

Introduction

Characterizing proteins and their post-translational modifications remain a major challenge for mass spectrometry. Electron capture dissociation (ECD) produces clean spectra of entire proteins while also preserving PTMs but remains impractical because of the technical difficulties of confining low-energy electrons. We have developed an efficient ECD technology that can be retrofitted into most common mass spectrometers.

Methods

The electromagnetostatic ECD cell is 30 mm long and consists of two magnets, an electron-emitting hot filament and 6 electrostatic lenses. The magnets provide a magnetic field that forces electrons emitted from the filament to be in the flight path of protein ions. The ECD cell is positioned between the mass resolving quad and the mass analyzers in either Orbitrap or in Q-ToF.

Results

We have developed an ECD cell that can be installed in one hour and does not affect ion transmission [1]. The cell was first optimized to maximize fragmentation of substance P. These same parameters worked well for fragmenting proteins up to 30 kDa. Sequence coverage exceeds 90% for carbonic anhydrase (29 kDa). In addition to producing c and z fragments, a substantial number of d and w sidechain fragments were present in the same spectrum that could distinguish many isobaric leucines from isoleucines. Labile post-translational modifications are retained. For example, all eight phosphoserines in alpha-casein (23 kDa) were retained with 87% sequence coverage. With superoxide dismutase, copper and zinc cofactors remained bound to c and z ions containing their respective binding sites. Deuterium labeling of ubiquitin enabled top-down hydrogen/deuterium exchange with rates consistent with prior NMR results. Sequence coverage of 80-95% for proteins like ubiquitin and a-synuclein (14 kDa) was obtained during UHPLC separations from peaks lasting 1-3 seconds.

Conclusions

The ECD cell is easily adapted into both Orbitrap and Q-ToF instruments without reducing their sensitivity or transmission. Fragmentation produces clean spectra that allows nearly complete characterization of proteoforms with preservation of labile PTMs. The presence of numerous d and w ions enhances the success of de novo sequencing by allowing many leucine/isoleucine or lysine/glutamine pairs to be easily distinguished. The ECD cell can help map disulfides in bottom up approaches [2].

Novel Aspect

High sequence coverage of proteins up to at least 30 kDa can be produced at chromatographic speeds by retrofitting the ECD into current mass spectrometers.

References

1. Voinov VG, Deinzer ML, Beckman JS, Barofsky DF. *J. Am. Soc. Mass Spectrom.* 22, 607-11 (2011)
2. Fort KL, Cramer CN, Voinov VG, Vasil'ev YV, Lopez NI, Beckman JS, Heck AJR. *J. Proteome Res.* 17, 926-933, (2018)

TP-170 / THE ROLE OF H₂O IN THE DISSOCIATION CHEMISTRY OF [W₆O₁₉]²⁻: A COMBINATIONAL STUDY BY MASS SPECTROMETRY AND DFT CALCULATION

Jie Cao⁽¹⁾ - *Chenggang Ci*⁽²⁾ - *Chang Liu*⁽¹⁾ - *JieChao Zang*⁽²⁾

Laboratory of Cluster Science, Ministry of Education of China, Beijing Key Laboratory of Photoelectronic/Electrophotonic Conversion Materials, School of Chemistry and Chemical Engineering, Beijing Institute of Technology, Beijing 100081, P. R. China (1) - School of Chemical and Chemical Engineering, Key Laboratory of Computational Catalytic Chemistry of Guizhou Province, Qiannan Normal University for Nationalities, Duyun 558000, P. R. China (2)

Keywords: Polyoxometalate; Tandem mass spectrometry; Dissociation mechanism; DFT calculation

Introduction

Polyoxometalates (POMs) are discrete anionic metal-oxygen clusters which exhibit a great diversity of applications in catalysis, medicine, etc. [1] The fragmentation mechanism of discrete polyanions as well as its protonated/organic ammonium/alkali metal adduct ions are well understood by the contribution of a few scientists [2-4], in contrast, the solvated polyanions (in particular H₂O as a solvent) chemistry is scarcely investigated. This sort of research can allow us to shed light on the mechanism of the inverse process of nucleation reaction in a solution, which attracts the most attention of synthetic chemists of POMs. Herein, we report our preliminary result regarding the water-solvated polyanions by taking the Lindqvist-type [W₆O₁₉]²⁻ as an example to probe the dissociation chemistry of the water-solvated [W₆O₁₉-H₂O]²⁻ in comparison to the bare [W₆O₁₉]²⁻ in order to explicitly define the role of water in the dissociation chemistry of polyanions by tandem mass spectrometry combined with density functional theory (DFT) calculation. Tandem mass spectrometry experiments of [W₆O₁₉]²⁻ and [W₆O₁₉-H₂O]²⁻ were performed on Agilent 6520 Q-TOF mass spectrometer with N₂ as the collision gas. The MS/MS results show that these two ions went through completely different dissociation pathways, where the former produces [W₃O₁₀]²⁻ and [W₄O₁₃]²⁻ fragments whilst the later generates predominately the mono-nuclear fragment [HWO₄]⁻. This result prompts us to search further about the reaction pathways for these two precursor ions by DFT calculation. The DFT results show that the reaction may proceed via two distinctive pathways: "water-free" and "water-assisted". The water-free pathway starts with the cleavage of W-Ob bond of [W₆O₁₉]²⁻ without transition state. Water-free process is calculated to be endothermic by about 39.1 kcal mol⁻¹. However, initially, in the water-assisted process, explicit water coordinates with POM by the intermolecular hydrogen bond, which reduces the energy for about 12 kcal mol⁻¹. Then the cleavage of W-Ob bond in the water-assisted pathway is found to be endothermic by about 30.1 kcal mol⁻¹. Next, in the water-free process, the cleavage of W-O bond is favorable thermodynamically, forming [WO₄]²⁻ anion and intermediate W₅O₁₅. However, for the water-assisted pathway, the water addition proceeds with an 8.2 kcal mol⁻¹ energy barrier at a transition state TS1. All calculated bond parameters and IRC calculation results indicate that this is a synchronous process, that is, the addition of H₂O accompanied by a H-transfer to form intermediate [HWO₄-W₅O₁₅OH]²⁻. Next the formation of [HWO₄]⁻ and [W₅O₁₅OH]⁻ proceeds without an energy barrier and is calculated to be exothermic by about 80.1 kcal mol⁻¹. The calculation is on the way to find out the lowest-energy routes to reach the final products, respectively.

Methods

Tandem mass spectrometry combined with density functional theory (DFT) calculation.

Results

Tandem mass spectrometry experiments of [W₆O₁₉]²⁻ and [W₆O₁₉-H₂O]²⁻ were performed on Agilent 6520 Q-TOF mass spectrometer with N₂ as the collision gas. The MS/MS results show that these two ions went through completely different dissociation pathways, where the former produces [W₃O₁₀]²⁻ and [W₄O₁₃]²⁻ fragments whilst the later generates predominately the mono-nuclear fragment [HWO₄]⁻. This result prompts us to search further about the reaction pathways for these two precursor ions by DFT calculation.

The DFT results show that the reaction may proceed via two distinctive pathways: “water-free” and “water-assisted”. The water-free and water-assisted pathway all start with the cleavage of W-Ob bond of $[W_6O_{19}]^{2-}$ (1) without transition state. Water-free process is calculated to be endothermic by about 39.1 kcal mol⁻¹. Similarly, the water-assisted pathway is found to be endothermic by about 30.1 kcal mol⁻¹. Next, in the water-free process, the cleavage of W-O bond is favorable thermodynamically, forming $[WO_4]^{2-}$ anion and intermediate W_5O_{15} . However, for the water-assisted pathway, the water addition proceeds with an 8.2 kcal mol⁻¹ energy barrier at a transition state TS1. All calculated bond parameters and IRC calculation results indicate that this is a synchronous process, that is, the addition of H₂O accompanied by a H-transfer to form intermediate $[HWO_4-W_5O_{15}OH]^{2-}$. Next the formation of $[HWO_4]^-$ and $[W_5O_{15}OH]^-$ proceeds without an energy barrier and is calculated to be exothermic by about 80.1 kcal mol⁻¹. The calculation is on the way to find out the lowest-energy routes to reach the final products, respectively.

Conclusions

The Lindqvist-type bare polyanion $[W_6O_{19}]^{2-}$ goes through completely different dissociation pathway from its water-solvated counterpart $[W_6O_{19} \cdot H_2O]^{2-}$, where $[W_6O_{19}]^{2-}$ produces $[W_3O_{10}]^{2-}$ and $[W_4O_{13}]^{2-}$ fragments whilst $[W_6O_{19} \cdot H_2O]^{2-}$ generates predominately the mono-nuclear fragment $[HWO_4]^-$. The preliminary result of DFT calculation shows that H₂O can largely lower the activation energy barrier of the initial steps of fragmentation.

Novel Aspect

This is the first investigation of the bare and the solvated polyanions with the objective to explicitly define the role of water in the dissociation chemistry of POMs.

References

1. M. T. Pope, *Heteropoly and Isopoly Oxometalates*; Springer-Verlag: New York, 1983
2. *Polyoxometalates: from Platonic Solids to Anti-Retroviral Activity*; M. T. Pope, A. Müller, Eds. ; Kluwer Academic Publishers: Dordrecht, The Netherlands, 1994.
3. Ma, M. T., Waters, T., Beyer, K., Palamarczuk, R., Richardt, P. J. S., O’Hair, R. A. J., Wedd, A. G., *Inorg. Chem.* 48, 598-606 (2009).
4. Gunaratne, K. D. D., Prabhakaran, V., Johnson, G. E., Laskin, J., *J. Am. Soc. Mass Spectrom.* 26, 1027-1035 (2015)
5. Johnson, G. E., Al Hasan, N. M., Laskin, J., *Int. J. Mass Spectrom.* 354-355, 333–341 (2013).
6. Cao, J., Li, C. C., Zhang, Z. X., Xu, C., Yan, J., Cui, F. Y., Hu, C. W., *J. Am. Soc. Mass Spectrom.* 23, 366-374 (2012)
7. J. Cao, C. Xu, Y. X. Fan, L. Y. Fan, X. H. Zhang, C. W. Hu, *J. Am. Soc. Mass Spectrom.*, 24, 884-894(2013).

TP-171 / CONTROLLED COLLISION INDUCED DISSOCIATION OF LARGE AROMATIC IONS IN A LQIT

Mickaël Carlos - Hassan Sabbah - Christine Joblin

Institut de Recherche en Astrophysique et Planétologie, Université de Toulouse, CNRS, CNES, TOULOUSE

Keywords: linear ion trap, PAH, collisional dissociation, ion dynamics, astrochemistry

Introduction

The AROMA (Astrochemistry Research of Organics with Molecular Analysis) setup [1] has been designed to analyse the organic molecular content of cosmic dust analogues, as part of the NANOCOSMOS ERC project. In this work, we test the capabilities of AROMA to study the energetics of molecular ion fragmentation using a Collision Induced Dissociation (CID) scheme. A motivation is to discriminate structural isomers.

Methods

AROMA is a hybrid setup, composed of a laser desorption/ionization source, coupled with a Linear Quadrupole Ion Trap (LQIT) and an orthogonal Time Of Flight (oTOF). Controlled CID experiments are performed in the LQIT under various conditions of colliding gas pressure and Dipolar Excitation (DE) conditions. Ion dynamics simulations using the Simion software are performed and combined with a Monte Carlo kinetics program.

Results

We present results on the CID of the pyrene cation $C_{16}H_{10}^+$ colliding with Ar gas. We explored the influence of the DE parameters (amplitude, frequency, duration of the DE). This includes the two different modes of excitation (on and off resonance). With these parameters we can tune the collisional heating rate, reaching different regimes in the competition between heating and cooling by dissociation. This competition governs the breakdown curve (loss of the parent cation) and impacts the H over 2H branching ratio.

Using Simion, we performed simulations of the ion motion taking into account collisional cooling by He during transfer to the LQIT and subsequent dipolar excitation in the LQIT. A hard sphere model was used to describe collisions. We obtained the distributions of the energies at the centre of mass (ECM) and of the time between collisions. Both distributions are used as inputs in our Monte Carlo model. Comparison between simulations and experiments allows us to constrain the efficiency of kinetic to internal energy conversion.

Conclusions

We have set a protocol to quantify fragmentation in the lowest energy channels. This can be used to discriminate structural isomers in complex mixtures such as cosmic dust analogues or extraterrestrial materials. In particular we have been able to reevaluate the assignment of pyrene as the carrier of $m/z=202.08$ in the Murchison meteorite [1]. In the future, we will demonstrate how our heating scheme can also be applied to more strongly bonded systems such as C_{60}^+ .

Novel Aspect

Our combined experimental and simulation approach allows us to better control dissociation at threshold for large aromatic ions submitted to collisions with Ar in a LQIT.

References

1. Sabbah, H., Bonnamy, A., Papanastasiou, D., Cernicharo, J., Martín-Gago, J. -A., Joblin, C., *Astrophys. J.*, 843, 34 (2017).

TP-172 / ULTRAVIOLET PHOTODISSOCIATION OF PEPTIDES AND PROTEINS, GENERATED BY ELECTROSPRAY AND MATRIX-ASSISTED LASER DESORPTION/IONIZATION, ON THE SAME Q-EXACTIVE MASS SPECTROMETER

Marialaura Dilillo⁽¹⁾ - **Erik de Graaf**⁽¹⁾ - **Avinash Yadav**⁽¹⁾ - **Mikhail Belov**⁽²⁾ - **Liam McDonnell**⁽¹⁾

Fondazione Pisana per la Scienza - ONLUS, Proteomics Laboratory, San Giuliano Terme (1) - Spectroglyph LLC, Spectroglyph LLC, Kennewick (2)

Keywords: UVPD, MALDI – UVPD, EP-MALDI – Orbitrap, Top-Down proteomics

Introduction

Recent literature reports demonstrate that ultraviolet photodissociation can provide extensive sequence coverage for ESI generated protein ions. UVPD has a much lower charge state dependence than ETD and CID, thus the singly charged protein ions produced by MALDI may exhibit equal propensity to fragment. Here we compare UVPD of MALDI and ESI generated ions using a Q-Exactive equipped with a EP-MALDI/ESI dual source and UVPD.

Methods

All experiments were performed using a Q-Exactive Plus (Thermo Scientific, Bremen) equipped with a ESI/EP-MALDI dual ion source (Spectrograph LLC, Kennewick WA, USA). UVPD was performed using a 193 nm ArF laser (Excistar XS 500, Coherent, Santa Clara CA, USA) aligned coaxially with the HCD cell.

ESI was performed using 10 μ M standard protein solutions infused at 5 μ L/min. MALDI samples were prepared using the dried droplet method (CHCA matrix, 1000:1 molar excess).

Results

UVPD of MALDI generated ions was benchmarked against MALDI-MSD and ESI-UVPD performed on the same accurate mass, high resolution mass spectrometer.

ESI was first used with peptide and protein standards to optimize laser trigger timing, number of laser pulses and laser pulse energy. UVPD of ESI generated ions using the Q-Exactive Plus performed identically to those previously reported using hybrid instruments.

The optimized parameters were used for UVPD of MALDI generated ions, specifically the small peptide bradykinin and the small protein thymosin β 4 (monoisotopic masses 1059.57 and 5205.63, respectively). The UVPD fragments were compared at different laser pulse energies, and with those obtained from HCD and MS/MS.

Conclusions

The results demonstrate that MALDI UVPD outperformed MALDI HCD for all laser pulse energies. While MALDI MS/MS also provided extensive sequence coverage, its lack of precursor ion selection limits its application to isolated proteins whereas MALDI UVPD can be used on proteins from complex mixtures.

Novel Aspect

UVPD performed on intact protein ions generated by EP-MALDI on a modified Q Exactive Plus, demonstrating its ability to further improve MS/MS capabilities for MALDI generated ions.

TP-173 / GAS PHASE REACTIONS OF INDOCYANINE GREEN USING FEMTOSECOND-LASER-PULSE INDUCED PHOTODISSOCIATION AND COLLISION-INDUCED DISSOCIATION

Elena Mitrofanov - Tassilo Muskat - Jürgen Grotemeyer

Institute of Physical Chemistry, University of Kiel, Kiel

Keywords: Photodissociation, Indocyanine Green, collision-induced dissociation, NIR

Introduction

The dye Indocyanine Green (ICG) is widely used for medical purposes and as a fluorescent dye [1]. Investigations on the fluorescence of ICG showed very short fluorescence lifetimes [2] which leads to the requirement of femtosecond laser pulses for photodissociation experiments. The results will be compared to collision-induced dissociation processes to investigate stability and reactions of ICG.

Methods

The mass spectra have been recorded on an Apex III FT-ICR mass spectrometer (7.05 T) (Bruker Daltonik, Germany) using different fragmentation techniques. The ionization was performed via Electrospray Ionization from different solvent mixtures. For the fragmentation reactions, collision gas (Argon) or a femtosecond laser system (Quantronix, USA) with 790 nm and 150 fs pulse duration were used.

Results

Depending on the solvent mixture, the sodiated or the protonated ions and their clusters were observed in the mass spectra. Collision-induced dissociation (CID) and photodissociation (PD) were applied to the different ions. The photodissociation process at a wavelength of 790 nm resulted in several intense

fragments for the sodium adduct and less intense fragments for the protonated molecule. Most fragments arise from the dissociation of a linkage between two identical aromatic systems. Beside these cleavages, only the bisodiated ICG shows the elimination of some small neutral fragments. In contrast to the PD spectra the CID-spectra showed poor fragmentation for the protonated molecule and no fragmentation for the sodium adduct.

Conclusions

The poor fragmentation in the CID process leads to the conclusion that the collision energy does not exceed the required energy that is reached by multiphoton absorption of the femtosecond laser pulses. The comparison of the fragments from the protonated and the sodiated molecules gives a hint for the influence of the sodium cation on the molecular stability. [3]

Novel Aspect

Photodissociation reactions of ICG and its clusters in the gas phase using femtosecond laser irradiation.

References

1. Alander J. T., Kaartinen I., Laakso A. et al., Int. J. Biomed. Imaging, 2012, ID 940585 (2012). DOI:10.1155/2012/940585.
2. Lee H., Berezin M. Y., Henary M. et al., J. Photochem. Photobiol. A, 200, 438 (2008). DOI: 10.1016/j.jphotochem.2008.09.008.
3. Mitrofanov E., Muskat T., Grotemeyer J., EJMS, published online 2017, in press 2018. DOI: 10.1177/1469066717750044.

TP-174 / LOST IN PHASE SPACE: FRAGMENTATION ENERGETICS IN QUADRUPOLE ION TRAPS

Thomas Neugebauer - Thomas Drewello

Friedrich-Alexander-University Erlangen-Nuremberg, Chair of Physical Chemistry I, Erlangen

Keywords: CID Fragmentation energetics, Quadrupole ion trap, Phase-space ellipses

Introduction

Quadrupole ion traps play a central role in the investigation of ion-fragmentation processes due to their capability to conduct MS_n experiments. The motion of ions inside quadrupole ion traps is complex but can be describes by Mathieu's equation. [1]Due to the complexity of the consequences of the change in stability parameters the fragmentation has been investigated in the past only by simulation programs. We present a systematical investigation of this process.

Methods

The analytical solution of Mathieu's equation has been combined with phase-space analysis. Phase-space ellipses have been implemented in calculations over their ellipse angle and explicitly not the quadratic equation, this way avoiding the common treatment of the ellipses as uniformly populated, thus providing a correct treatment of time.

Results

The phase-shift between the fundamental ion oscillation and the superimposed micromotion creates an elliptic relationship between the position and the velocity of ions at a specific phase of the applied alternating potential. We can show by appropriate application of trigonometry that the angle (the position on the ellipses in phase-space has regarding the coordinate system axes)is related to the variable time and with that, the probability to find an ion on a certain position in phase-space varies.

By the use of these ellipses we provide a straightforward explanation, why and how the amplitudes of ion

trajectories can drastically change, resulting in loss of ions, although their m/z value is above the Low-mass cut-off (LMCO) and to which extent this is expected. [2] By the use of the center-of-mass collision energy, E_{com} , we define effective temperatures and further give information about the likelihood of ions enduring successive excitation that can result in consecutive fragmentation. [3]

Conclusions

We explain how the variables are connected to each other in the complex fragmentation process and provide clues regarding the origin of certain observations as being caused by instrumental or molecular properties. Further general trends are given and it is also discussed how a change in LMCO settings might help achieving the desired information.

Novel Aspect

Consequences of the change in stability parameters during fragmentation are presented with special emphasis on the likelihood for successive fragmentation events.

References

1. Paul W., Ang. Chem. Int. Ed., 29, 739-748 (1990)
2. T. S. Neugebauer, T. Drewello (in preparation)
3. T. S. Neugebauer, T. Drewello (in preparation)

TP-175 / COLLISION INDUCED DISSOCIATION OF HYDROGEN DEFICIENT RADICAL IONS ON THE OMNITRAP PLATFORM

Dimitrios Papanastasiou⁽¹⁾ - Roman A. Zubarev⁽²⁾

Fasmatech Science & Technology, Patriarchou Gregoriou, Agia Paraskevi, Athens (1) - Karolinska Institutet, Department of Medical Biochemistry and Biophysics, Stockholm (2)

Keywords: top-down, ion activation, linear ion trap

Introduction

The omnitrapp platform is a novel linear ion trap sustaining a diverse ion activation network including collision-induced and electron directed dissociation of ions. A series of experiments involving collision induced dissociation of radical ions formed by meta-ionization of poly-protonated species using $>10\text{eV}$ electrons are performed and new informative transitions in the fragmentation patterns are reported for insulin chain B and ubiquitin ions.

Methods

Precursor ions selected in the quadrupole of a QExactive+ are transferred into the omnitrapp platform. Electron meta-ionization of protonated ions is performed by varying electron energy, adjusting irradiation time and the frequency of the trapping waveform to maximize radical ion abundance. CID MS3 experiments of EI-produced hydrogen-deficient radicals are compared to their protonated counterparts in the same charge state.

Results

Electron meta-ionization of insulin chain B ions is performed and subsequent CID of the radical ions reveals the formation of three complementary c, z fragment pairs near tyrosine residues. A mechanism is proposed to explain the dissociation pathway involving meta-ionization at or near the tyrosine residue followed by the rapid phenolic proton migration towards a basic site, *e.g.* N-terminus, forming a protonated group, with a neutral radical remaining at the tyrosine residue and initiating N-Ca backbone cleavage. Electron meta-ionization of 7^+ ubiquitin ions produced $8^+(\cdot)$ hydrogen deficient species. Comparison of CID MS/MS data of the 8^+ protonated ions with the hydrogen deficient radicals at the same charge state reveal significant

differences, especially by the presence and position of a, x and c, z fragmentation patterns, as well as significant similarities in the relative abundances of b, y fragments.

Conclusions

The unique versatility of the ion activation network available in the omnitrapp platform is demonstrated using MS³ CID experiments of meta-ionized species. Enhanced sequence information is obtained in CID experiments of the radical ions compared to even electron species at the same charge state. Cleavage position, fragment type and abundance are analyzed in great detail revealing the hitherto unreported fragmentation processes.

Novel Aspect

Novel ion activation tools for the orbitrap mass analyzer and expanded versatility of the omnitrapp platform for top-down mass spectrometry.

TP-176 / COMPREHENSIVE GLYCOMICS STUDY BY NANOLC-EED-MS/MS

Yang Tang^(1,2) - ***Juan Wei***⁽²⁾ - ***Catherine Costello***^(1,2) - ***Cheng Lin***⁽²⁾

Department of Chemistry, Boston University, Boston, MA (USA) (1) - Department of Biochemistry and Center for Biomedical Mass Spectrometry, Boston University School of Medicine, Boston, MA USA (2)

Keywords: Glycomics, PGC-LC, EED-MS/MS

Introduction

Comprehensive structural elucidation of glycan mixtures is challenging. Here, glycans are reduced and permethylated or labeled with a reducing-end fixed charge to increase sensitivity, avoid gas-phase structural rearrangement and facilitate spectral interpretation. A PGC column was used to maximize isomer separation [1]. EED was utilized due to its ability to generate informative MS/MS spectra to allow accurate determination of both the topology and linkages.

Methods

Glycan standards, LNT, LNnT, LNFP I, II, III, V and VI, were either reduced and permethylated [1], or labeled with a reducing-end fixed charge via reductive amination with 3-aminopyridine followed by methylation using iodomethane in ACN [2]. Separation was carried out on a Waters nanoACQUITY UPLC, with a nanoPGC column held at 50 °C. EED-MS/MS analyses were performed on a Bruker 12-T solarix FTICR MS with the cathode bias set between 16 V and 18 V [3-4].

Results

PGC-LC can often achieve better isomer resolution than RPLC. For example, reduced and permethylated linkage isomers, LNFP II and III, are barely separated by RPLC, yet they are baseline resolved by PGC. Meanwhile, glycans with a reducing-end fixed charge are not retained by RPLC, but form sharp peaks in PGC.

EED can generate more informative spectra than CID or HCD. Here, for glycans with a reducing-end fixed charge, EED produced complete series of Z-, Y- and 1,5X-ions, forming triplets with characteristic spacing that could be easily recognized to facilitate accurate determinations of the glycan topology. Additional radical-driven dissociation pathways were also identified, generating linkage-specific cross-ring, secondary or internal fragments throughout the glycan backbone. The presence and location of the branching site(s) and the linkage configuration at HexNAc residues could also be determined based on characteristic fragmentation patterns occurring only at those sites.

Conclusions

Isomer resolution was achieved on PGC-LC for glycans in two different derivatized forms. EED of

permethylated glycans and glycans with a reducing-end fixed charge produced rich structural information with unique spectral features for reliable determination of the glycan sequence, branching pattern and linkage configurations. Online PGC-LC-EED-MS/MS should become a powerful method for comprehensive glycome characterization.

Novel Aspect

Separation and detailed structural characterization of isomeric glycan mixtures are achieved in a single LC-EED-MS/MS analysis.

References

1. Costello C. E., Contado-Miller J. M., Cipollo J. F., *J. Am. Soc. Mass. Spectrom.*, 18, 1799-1812 (2007)
2. Hase S., Ikenaka T., Matsushima Y., *Biochem. Biophys. Res. Commun.*, 85, 257-263 (1978)
3. Yu X., Huang Y., Lin C., Costello C. E., *Anal. Chem.*, 84, 7487-7494 (2012)
4. Yu X., Jiang Y., Cheng Y., Huang Y., Costello C. E., Lin C., *Anal. Chem.*, 85, 10017-10021 (2013)

TP-177 / ESI(-)MS/MS OF 2-OXOCARBOXYLATE HYDRAZONES: CARBOXYL GROUP PARTICIPATION IN TWO FRAGMENTATION PATHWAYS

Robert L. White - Ben W. Firth - Luc M. LeBlanc - J. Stuart Grossert

Dalhousie University, Chemistry, Halifax

Keywords: ESI(-)MS/MS, hydrazone derivatives, fragmentation mechanism, decarboxylative elimination, computations

Introduction

Analysis of aldehydes and ketones has been achieved by formation of hydrazone derivatives and negative-ion electrospray tandem mass spectrometry [ESI(-)MS/MS] [1,2]. Upon collision induced dissociation (CID), major product ions were formed by losses of NO and NO₂ [3]. In our current ESI(-)MS/MS studies on hydrazones of 2-oxo carboxylic acids, fragmentations are directed by the carboxylate group.

Methods

Mass spectra were obtained using ESI(-)MS/MS with MeOH solvent. A triple quadrupole spectrometer gave both product- and precursor-ion spectra, while an ion trap spectrometer gave complementary product-ion spectra. Quantum mechanical computations of gas-phase ions were done using Gaussian 09 software at the MP2/6-311++G(2d,p)//B3LYP6-31++G(2d,p) level of theory.

Results

Phenyl, 4-nitrophenyl and 2,4-dinitrophenylhydrazone derivatives of glyoxylic and phenylglyoxylic acid were prepared. Upon CID, the deprotonated glyoxylate and phenylglyoxylate hydrazones primarily yielded anilide ions; decarboxylation was a minor process for the glyoxylate derivatives. The major pathway was consistent with decarboxylative elimination [4], in which an anilide ion, a nitrile and CO₂ are formed by N–N and C–C bond cleavage. Barriers of less than 100 kJ/mol were determined by computations and breakdown curves showed the influence of nitro and phenyl substituents on the energetics of decarboxylation. ESI of samples in acetonitrile gave reduced abundance of the minor decarboxylated product ion. The minor pathway was attributed to solvent-assisted ionization on nitrogen, rearrangement to Ar-N=N-CH₂CO₂⁻ (barrier = 130 kJ/mol) and lastly decarboxylation (negligible barrier).

Conclusions

The fragmentation behavior of the 2-oxo carboxylate hydrazones resulted from reactions of the carboxyl group. The major fragmentation pathway proceeded by decarboxylative elimination, whereas the minor pathway observed for the glyoxylate hydrazones involved ionization on nitrogen, rearrangement and

subsequent decarboxylation.

Novel Aspect

Different pathways of deprotonated hydrazone derivatives highlight the challenges of predicting the fragmentation behavior of multifunctional ions.

References

1. Szulejko J. E., Kim K-H., Trends Anal. Chem., 64, 29 (2015).
2. Uchiyama S., Inaba Y., Kunugita N., J. Chromatogr. B, 879, 1282 (2011).
3. Zwiener C., Glauner T., Frimmel F. H., Anal. Bioanal. Chem., 372, 615 (2002).
4. LeBlanc L. M., Powers S. W., Grossert J. S., White R. L., Rapid Commun. Mass Spectrom., 30, 2133 (2016).

TP-178 / ASSESSING BINDING AFFINITY WITH LILBID-MS

Genia Hense - Phoebe Young - Nina Morgner

Institute of Physical and Theoretical Chemistry, Goethe Universität Frankfurt, Frankfurt am Main, Germany

Keywords: binding affinity, dissociation, LILBID, native mass spectrometry, ion source

Introduction

LILBID is an ion desorption method for native mass spectrometry in which droplets of aqueous solution are irradiated with pulsed IR laser light, leading to explosive expansion of the droplets and release of sample molecules as ions. LILBID-MS is primarily used to assess the stoichiometry of noncovalent biomolecular complexes. Here we outline a new method to quantitatively determine binding affinities with LILBID-MS and present preliminary results from this method.

Methods

In LILBID, harsher or softer laser energy transfer leads to different amounts of dissociation of noncovalently bound biomolecular complexes. In this new method, the sample is measured while purposely varying the energy transfer from the desorption laser to the droplet in the LILBID ion source. The amount of dissociation in the resulting spectra is then plotted against a measure of energy input, resulting in a dissociation curve analogous to a melting curve.

Results

The amount of energy transferred to the sample during desorption can be changed by varying instrument parameters such as laser power and droplet position within the laser beam. However the relationship between these parameters and the amount of energy transferred to the sample is not straightforward. We monitor energy input using the shape of the droplet's explosive expansion. Certain parameters of the shape of this expansion are closely correlated with both energy input and spectral results and can be used as robust measures of laser energy transfer.

During a binding affinity measurement, one of the instrumental parameters for energy transfer is purposely varied as droplets are emitted. For each droplet, an explosion image and a mass spectrum are recorded. The shapes of the explosive expansions are analyzed, and the amounts of dissociation in the spectra are calculated and plotted against the corresponding laser energy transfer, which we derive from the explosion shape. The resulting curve is a plot of dissociation with increasing energy input.

Conclusions

Dissociation curves for different short dsDNA reflect the binding strengths of the samples, demonstrating that this method allows conclusions about the binding strength of biomolecules in solution. Since data for each dissociation plot can be collected within minutes and with a fraction of the sample necessary for a titration or a melting curve with LILBID, this method represents an attractive alternative to titration-based

binding affinity experiments.

Novel Aspect

This method constitutes a new way to perform fast top-down analysis of binding affinities of biomolecular complexes.

TP-179 / TUNING MASS SPECTROMETRIC PROTEOMICS EXPERIMENTS: COLLISION ENERGY DEPENDENCE OF MASCOT SCORE

Ágnes Révész⁽¹⁾ - *Tibor Andras Rokob*⁽²⁾ - *Dany Jeanne Dit Fouque*⁽³⁾ - *Lilla Turiák*⁽¹⁾ - *Antony Memboeuf*⁽³⁾
- *Károly Vékey*⁽¹⁾ - *László Drahos*⁽¹⁾

Research Centre for Natural Sciences, Hungarian Academy of Sciences, MS Proteomics Research Group, Budapest, Hungary (1) - Research Centre for Natural Sciences, Hungarian Academy of Sciences, Theoretical Chemistry Research Group, Budapest, Hungary (2) - Université de Bretagne Occidentale, UMR CNRS 6521, CEMCA, Brest, France (3)

Keywords: proteomics, peptide fragmentation, collision energy, optimization, Mascot score

Introduction

Recently, MS-based proteomics has become a powerful tool in protein identification and characterization. In the bottom-up approach, the choice of collision energy of the peptide MS/MS experiment significantly affects the performance of the method. High-level objectives, e.g. highest identification score, have not been explored in this respect. The present work maps the energy dependence of this proteomics variable and determines the optimal collision energy.

Methods

Tryptic digests of complex standards (HeLa, E. Coli) were subject to nano-HPLC-MS/MS experiments at 21 different collision energies on a Bruker Maxis II ETD QTOF instrument. Spectra were matched against the SwissProt database using the Mascot search engine. Further analysis, including comparison of LC-MS/MS runs and fitting of the Mascot score vs. collision energy curves, were performed with an in-house developed software using Mascot output files as inputs.

Results

A cluster analysis of Mascot score vs. energy curves revealed two types of energy dependent behavior (unimodal and bimodal) in about 50–50% of the investigated peptides. For bimodal peptides, number and summed intensity of b fragment ions show maxima around the first score maximum. y ions are less sensitive to collision energy but more intense and numerous near the second score maximum. The optimal collision energy, obtained as the peak positions from one or two fitted Gaussians, follow linear trends with respect to m/z , but the three groups have different m/z dependence. The presence of mobile proton decreases the optimum significantly (1–5 eV), but other factors also seem important indicated by the large residual variance even of the separate linear fits for peptides with mobile and partially mobile protons. Where 2 maxima exist, their difference is surprisingly m/z independent and it is ~ 17 eV on average for both +2 and +3 peptides.

Conclusions

An unprecedented bimodal behavior of the score vs. energy curves, linked to the different behavior of y- and b-type ions, was identified for about half of the peptides. For these, the average score increases by $\sim 15\%$ upon energy optimization. On the other hand, for unimodal peptides, the original default collision energy seemed quite good, only a slight increase of $\sim 4\%$ could be achieved. The optimal energy has linear m/z dependence with marked effect of proton mobility.

Novel Aspect

This is the first study to address collision energy dependence of a peptide identification score and to quantify the role of less abundant b-ions in peptide identification [1].

References

1. Revesz A., Rokob T. A., Jeanne DitFouque D., Turiak L., Memboeuf A., Vekey K., Drahos L., J. Proteome Res., 17, 1898-1906 (2018).

TP-180 / STRUCTURAL CHARACTERIZATION OF COMPLEX LIPIDS BY OZONE-INDUCED DISSOCIATION AND ULTRAVIOLET PHOTODISSOCIATION ON HIGH-RESOLUTION MASS SPECTROMETERS

Angela Criscuolo⁽¹⁾ - *David L. Marshall*⁽²⁾ - *Martin Zeller*⁽¹⁾ - *Berwyck L. J. Poad*⁽²⁾ - *Jan-Peter Hauschild*⁽¹⁾ - *Todd W. Mitchell*⁽³⁾ - *Gavin Reid*⁽⁴⁾ - *Stephen Blanksby*⁽²⁾

Thermo Fisher Scientific (Bremen) GmbH, Life Science Mass Spectrometry, Bremen, Germany (1) - Queensland University of Technology, Central Analytical Research Facility, Institute for Future Environments, Brisbane, Australia (2) - University of Wollongong, School of Medicine, Illawarra Health and Medical Research Institute, Wollongong, Australia (3) - The University of Melbourne, School of Chemistry, Department of Biochemistry and Molecular Biology, Parkville, Australia (4)

Keywords: Lipidomics, ultraviolet photon dissociation, ozonolysis, double-bond position

Introduction

Lipids play different physiological roles as sources of stored metabolic energy or structural matrices in biological membranes. Alterations to the lipidome have been associated with the onset and progression of numerous diseases such as cancer, diabetes and cardiovascular diseases. Detailed structural characterization of lipids in complex extracts is thus desirable, including the determination of: (1) lipid class; (2) number of carbons and double bonds; (3) relative position of acyl chains; and (4) location and stereochemistry of double-bonds.

Current tandem mass spectrometry approaches identify lipid class, total carbons and double bonds, but not further information. In this study, we investigate the utility of new ion activation techniques, specifically, ozone-induced dissociation (OzID) and ultraviolet photodissociation (UVPD) for more comprehensive lipid structure elucidations.

Methods

Lipids were extracted from human blood plasma [1] and analyzed using both shotgun and LC-MS approaches. For the latter, lipids were separated using a Thermo Scientific™ Accucore™ C30 column (2.1 x 150 mm, 2.6 μm) and a Thermo Scientific™ UltiMate 3000™ UHPLC system operated at 50°C, and a flow rate of 325 μL/min. ESI-MS analysis was performed in positive ion mode on a modified Q Exactive HF mass spectrometer with ozone (approximately 11% O₃ in oxygen O₂) replacing nitrogen as the HCD collision gas. Photodissociation was implemented on a Thermo Scientific™ Orbitrap Fusion™ Lumos™ Tribrid™ MS equipped with a 213 nm UVPD laser.

Results

Preliminary experiments were performed incorporating OzID on a Thermo Scientific™ Q Exactive™ HF MS to identify the position of acyl chain position and double-bonds in phospholipids and triacylglycerols. By combining HCD activation with OzID on selected [M + Na]⁺ lipid ions on the Thermo Scientific™ Q Exactive™ HF MS, the relative position of acyl chains on the glycerol backbone could be determined without compromising in scan time. In a shotgun phospholipidome analysis of plasma it was possible to cover the mass range from *m/z* 650 to 950 in only 14 minutes. Importantly, a single phospholipid HCD/OzID spectrum provides information about the phospholipid class, the number of carbons and double bonds and their sn-position, therefore highlighting the presence of isomeric lipids.

Additionally, reaction of selected lipid ions with ozone inside the mass spectrometer can produce two characteristic product ions: an aldehyde and a Criegee ion. These ions are diagnostic for the double bond position and can be obtained by OzID increasing the HCD cell fill time and decreasing the collision energy. With this technique, double bond position of the analyzed lipids can be readily obtained. We found the triacylglycerol and phospholipid fractions within plasma samples to be rich in monounsaturated fatty acids containing a double bond in either the n-7 or n-9 position. OzID experiments were compared with lipid fragmentation by UVPD on an Orbitrap Fusion Lumos MS equipped with a 213 nm laser. Experiments conducted by Williams et al. [2] and Ryan et al. [3] report this new tool as a promising technique to complete structural elucidation of lipid structures.

For Research Use Only.

Conclusions

In our study we compared ozonolysis and ultraviolet photon-dissociation for comprehensive lipid structure elucidations. Insaturation in cholesterol ester, phospholipids and triacylglycerols can be identified with both techniques due to the generation of structural diagnostic product ions: the aldehyde and a Criegee ion for the Ozonolysis and a pair of ions formed by the cleaving in specific site of the double bond upon photon uptake. Linoleic acid (18:2n-6, n-9) was observed as the major polyunsaturated fatty acid in the lipidome with both techniques.

Novel Aspect

Successful implementation of OzID and UVPD as complementary techniques for comprehensive global lipid structure elucidation and resolving isomeric lipids.

References

1. Matyash V. et al. J. Lipid Res., 2008, 49, 1137–1146
2. Williams et al. JACS 2017, 139, 15681-15690.
3. Ryan et al. J. Am. Soc. Mass Spectrom. 2017, 28, 1406-1419.

TPS-S09 – Astrochemistry

Chairs: Paul Mayer, Simonetta Fornarini

TP-181 / PIRENEA 2: A NEW SETUP TO STUDY COSMIC DUST IN THE LABORATORY

Anthony Bonnamy⁽¹⁾ - **Pavol Jusko**⁽¹⁾ - **Ming-Chao Ji**⁽¹⁾ - **Patrick Moretto-Capelle**⁽²⁾ - **Sebastien Zamith**⁽²⁾ - **Jean-Marc L'Hermite**⁽²⁾ - **Michel Pellarin**⁽³⁾ - **Christine Joblin**⁽¹⁾

Institut de Recherche en Astrophysique et Planétologie, Université de Toulouse, CNRS, CNES, TOULOUSE (1) - Laboratoire Collisions Agrégats Réactivité, Université de Toulouse, CNRS, TOULOUSE (2) - Institut Lumière Matière, Université Lyon 1, CNRS, LYON (3)

Keywords: Laboratory astrophysics, clusters, molecular aggregates, cryogenic ion trap, FTICR mass spectrometry

Introduction

The ERC Synergy Nanocosmos project aims to understand the physical and chemical processes leading to the formation of cosmic dust. In the framework of Nanocosmos, the PIRENEA 2 setup, an upgrade of the PIRENEA setup[1], is currently under construction in our laboratory. It combines sources of cosmic dust analogues with cryogenic trapping to study the interaction of relevant species with photons and gas under interstellar conditions.

Methods

The species that are generated from either a molecular aggregation source [2] or a laser vaporization cluster source [3], can be mass-selected and stored in a cryogenic radio-frequency (RF) quadrupole trap to be thermalized at tunable temperatures (~4-300 K). They can then be transferred to a cryogenic ion cyclotron resonance (ICR) cell (~20 K) for high-resolution mass analysis. Interaction with photons and gas can be studied in the quadrupole trap and in the ICR cell.

Results

The different parts of the setup (sources, RF trap and ICR cell) are still in the commissioning phase. In this poster, we will report on the achieved progress. Since polycyclic aromatic hydrocarbons (PAHs) are considered as the smallest in size dust particles, we will show how derivatives of these PAHs, which can be of astrophysical interest, can be produced with our sources including PAH clusters [4] and organometallic complexes made of PAHs and heavy atoms (Fe, Si) [5].

On the other hand, first tests of the performances of the new Fourier Transform ICR mass spectrometer (FTICR-MS) were performed. In particular this setup is equipped with a cryogenic pre-amplifier that was installed in order to significantly increase the sensitivity of ion detection [6].

Conclusions

Ion sources and the FTICR-MS part of PIRENEA 2 were tested. Next step will be to commission the cryogenic RF ion trap and interface all the modules together to get PIRENEA 2 fully operational.

Novel Aspect

The new PIRENEA2 setup opens new perspectives for the study in the laboratory of large molecular species and clusters under interstellar conditions.

References

1. Joblin C., Pech C., Armengaud M., Frabel P., Boissel P., EAS Publ. Ser., 4, 73 (2002).
2. Braud I., Zamith S., L'Hermite J. -M., Rev. Sci. Instrum., 88, 43102 (2017).
3. Truong N. X., Haertelt M., Jaeger B. K. A., Gewinner S., Schöllkopf W., Fielicke A., Dopfer O., Int. J. Mass Spectrom., 395, 1 (2016)
4. Rapacioli M., Calvo F., Joblin C., Parneix P., Toubanc D., Spiegelman F., Astron. Astrophys., 460, 519 (2006).
5. Simon A., Joblin C., J. Phys. Chem. A, 113, 4878 (2009).
6. Mathur R., Knepper R. W., O'Connor P. B., J. Am. Soc. Mass Spectrom., 18, 2233 (2007).

TP-182 / MASS SPECTROMETRY A NEW APPROACH TO STUDY THUNDERSTORMS

Luan Felipe C. Oliveira - David Ulisses Tega - Fernando Galembeck - Marcos N. Eberlin

UNICAMP, Chemistry Institute, Campinas

Keywords: Mass Spectrometry, Thunderstorms, Separation of Charge, V-EASI

Introduction

Until now in the study of thunderstorms techniques to measure voltage, charge, current, electromagnetic field and visible radiation were employed to understand about their electric aspects [1-5]. Alternatively, herein we present a precursor work in the use of mass spectrometry to understand the molecular aspects involved in thunderstorms. Therefore, physical changes in water were studied simulating some changes in clouds.

Methods

A quadrupole-mass spectrometer (LCMS Shimadzu 2010) with Ventury Easy Ambient Sonic-Spray ionization (V-EASI) was used in positive and negative mode with a range of m/z 10-1000 Da to analyze Milli Q water during melting process and warming in the range of 25°C – 80°C. Vapor of Milli Q water also measured without any source of ionization during warming in the range of 25°C – 80°C and after its boiling point.

Results

During the melting process different rates of variation in ionic current (IC) were observed between 10-50%. In the study of the vapor were observed an increase of the IC directly proportional of the amount of water vapor, the IC was ten times higher in this process. In the study of water warming a rise of two times in the IC only between 25°C – 60°C was present in the positive mode, between 60°C – 80°C the IC was constant. On the other hand, in the negative a rise in the whole interval studied was observed and an increase of four times. In all cases, we noted that the species responsible for the rise of IC were water clusters. The rise of charge during the melting is in agreement with other works[6-8], while the charge increase during the range of temperature and liquid-vapor transition is new among the processes of formation of clouds and thunderstorms.

Conclusions

This work presents a new application of Mass Spectrometry, the study of thunderstorms, showing to be efficient to study processes of electrification already studied and also new ones. The most important aspect of this work is that the technique provides information about molecular aspects involved in thunderstorms such as the suggestion that the atmospheric discharge could be a discharge of water clusters.

Novel Aspect

The use of mass spectrometry to study process related to thunderstorms.

References

1. Siingh D., Gopalakrishnan V., Singh R. P., Kamra A. K., Singh S., Pant V., Singh R., Singh A. K. (2007) The atmospheric global electric circuit: An overview. *Atmospheric Research*, 84, 91-110.
2. Rycroft M. J., Harrison R. G., Nicoll K. A., Mareev E. A. (2008) An Overview of Earth's Global Electric Circuit and Atmospheric Conductivity. *Space Science Reviews*, 137, 83-105.
3. Pinto I. R. C. A., Pinto O. (2003) Cloud-to-ground lightning distribution in Brazil. *Journal of Atmospheric and Solar-Terrestrial Physics*, 65, 733-737.
4. Farias W. R. G., Pinto O., Naccarato K. P., Pinto I. R. C. A. (2009) Anomalous lightning activity over the Metropolitan Region of São Paulo due to urban effects. *Atmospheric Research*, 91, 485-490.
5. Chiu, C. -S., Klett J. D. (1976) Convective electrification of clouds. *Journal of Geophysical Research*, 81, 1111-1124.
6. Carignano M. A., Shepson P. B., Szleifer I. (2007) Ions at the ice/vapor interface. *Chemical Physics Letters*, 436, 99-103.
7. Dash J. G., Mason B. L., Wettlaufer J. S. (2001) Theory of charge and mass transfer in ice-ice collisions. *Journal of Geophysical Research: Atmospheres*, 106, 20395-20402.
8. Saunders C. (2008) Charge Separation Mechanisms in Clouds. *Space Science Reviews*, 137, 335-353.

TP-183 / REACTION DYNAMICS OF SUBSTITUTED PAH IONS

Paul Mayer - Brandi West - Lukas Lesniak - Jake Burner

University of Ottawa, Dept. of Chemistry and Biomolecular Sciences, Ottawa

Keywords: polycyclic aromatic hydrocarbon ions, astrochemistry, iPEPICO, collision-induced dissociation

Introduction

The chemical origins of the Diffuse Interstellar Bands have remained elusive. Polycyclic aromatic hydrocarbons, PAHs have been postulated to be sources of DIBs due to the similarity between the IR emissions of the cations of PAHs and the DIBs. One molecule that has been investigated for its potential role as a DIB source is 1-methylpyrene. In this study we examine the dissociation of this ion and two other substituted PAH ions, 1-nitropyrene and 9-cyanophenanthrene.

Methods

iPEPICO experiments were conducted on the VUV beamline at the Swiss Light Source (Paul Scherrer Institut, Villigen, Switzerland). The mass-selected TPES were acquired to establish the IEs of each compound. RRKM theory was used to determine the 0 K activation energies (E_0) and entropies of activation ($\Delta^\ddagger S_{1000K}$) by modeling the experimental breakdown diagrams. Reaction pathways were calculated at the UCCSD/6-31G(d)//B3-LYP/6-31G(d) level of theory.

Results

For ionized 1-methylpyrene, the loss of H could form either 1-methylenepyrene (A) or a tropylium-motif containing ion (B). The transition state of highest energy in the pathway to B lies ~ 3.20 eV, while that for A is only 2.54 eV, favoring simple H loss from the methyl group, in contrast to that observed for ionized toluene. The presence of the three additional rings in 1-methylpyrene prevents the 7-member ring from expanding to relax the structure and lower the TS energy.

For ionized 1-nitropyrene, competition was observed between NO loss (resulting from a nitro-nitroso isomerization in the precursor ion), the loss of NO_2 to form a triplet state product ion, and a sequential loss of CO, which was confirmed by D-labelling experiments.

Unlike the previous 2 ions, the $-\text{CN}$ group in 9-cyanophenanthrene appears to act only to destabilize the PAH product ions. The $-\text{CN}$ group itself is not directly involved in the dissociation, shown to be H and C_2H_2 loss. The energies required for these reactions (~ 5.5 eV) are higher than observed for unsubstituted PAHs (~ 4.5 eV, regardless of PAH size) [1], an effect traceable to the impact of the electron-withdrawing $-\text{CN}$ group on the product ion stability.

Conclusions

The substituent on a large PAH like pyrene has been shown to play a variety of roles in the unimolecular dissociation of the ions. It can either drive the reaction mechanism (nitropyrene) or have a predominantly electronic effect (cyanophenanthrene). The larger ring system of pericondensed PAHs can also push the chemistry away from extensive isomerization (methylpyrene).

Novel Aspect

Presented are the first detailed investigations of the unimolecular ion chemistry of 3 substituted PAHs.

References

1. West B, Rodriguez Castillo S, Sit A, Mohamad S, Lowe B, Joblin C, Bodi A, Mayer PM, Phys. Chem. Chem. Phys., 20, 7195 (2018).

TP-184 / A MULTI-VENDOR APPLICABLE TARGET SYSTEM FOR (MA)LDI INSTRUMENTS WITH DISTINCT PERFORMANCE CHARACTERISTICS FOR THE ANALYSIS OF SMALL ORGANIC MOLECULES RELATED TO C CHONDRITES

Edita Rados⁽¹⁾ - ***Ernst Pittenauer***⁽¹⁾ - ***Johannes Frank***⁽²⁾ - ***Kurt Varmuza***⁽³⁾ - ***Günter Allmaier***⁽¹⁾

TU Wien, Institute of Chemical Technologies and Analytics, Vienna (1) - TU Wien, Joint Workshop of the Faculty of Technical Chemistry, Vienna (2) - TU Wien, Institute of Statistics and Mathematical methods in Economics, Vienna (3)

Keywords: small organic molecules, LDI, target adapter, gold target

Introduction

A serious comparison of one and the same sample preparation on different LDI-MS instruments is usually not possible because each device requires the use of its own, specific target. Therefore, we developed a target system, which enables the use of only one target for different devices. The work was focused on the analysis of small organic molecules selected due to their different desorption/ionization behavior and their potential presence in carbonaceous chondrites.

Methods

The selected organic molecules (*e.g.* tryptophan, 2-deoxy-D-ribose and triphenylene) were prepared and spotted onto the target plate (stainless steel or gold) and analysed on three devices with different performance characteristics, namely an Axima TOF2 (LinTOF/curved field RTOF), an ultrafleXtreme (LinTOF/dual stage RTOF) and a Synapt G2 (Qq/dual stage RTOF). The samples were also analysed on the TOF-SIMS 5 mass spectrometer equipped with gridless RTOF.

Results

The developed target system consisting of modified Bruker target adapter, a standard SS Waters target plate and an Au target plate allowed comparison of data acquired on four MS instruments with different performance characteristics (*e.g.* different lasers, different laser repetition rates, different ion source vacuum regimes). The data were obtained for all analytes from plain SS and Au target. In positive-ion mode, we were able to detect different types of molecular ions (M^+ , radical cations, $[M+H]^+$, $[M+Na]^+$, and $[M+2Na-H]^+$ ions) depending on the type of analyte. All sample molecules showed similar desorption/ionization behavior on all instruments in the positive-ion mode independently of the target applied. A deprotonated molecule of tryptophan was also detectable on all instruments with the use of SS and Au target. The use of the Synapt G2 instrument allowed accurate mass determination ($<\pm 4$ ppm) at a resolution of up to RFWHM 20.000 for all types of molecular ion. Additionally, low- and high-energy CID-spectra were acquired for tryptophan.

Conclusions

We showed that an in-depth comparison of the data obtained with devices from different vendors with distinct performance characteristics is possible when the presented target system is being used. The possibility to obtain accurate masses can be useful in identifying unknown compounds in C chondrite meteorite samples.

Novel Aspect

The developed target system allows the use of only one target on multiple LDI/MALDI-devices. Small organic molecules can be measured with LDI and accurate mass can be obtained.

TPS-S10 - Life sciences - Metabolomics

Chairs: Ulf Sommer, Donatella Caruso

TP-185 / MEASURING TREHALOSE 6-PHOSPHATE (TRE6P) IN DIFFERENT PLANT MATERIALS – THE TECHNICAL CHALLENGE OF WORKING AT THE FEMTOMOLE SCALE

Regina Feil - Maria Grazia Annunziata - John E. Lunn

Max Planck Institute of Molecular Plant Physiology, Dept. Metabolic Networks, Potsdam-Golm, Germany

Keywords: Trehalose 6-phosphate (Tre6P), High-Performance Anion-Exchange Chromatography (HPAEC), Tandem mass spectrometry (MS/MS), Stable-isotope labelled internal standards (SIL-IS).

Introduction

Trehalose 6-phosphate (Tre6P) is an essential metabolite in plants. It acts as a signal of sucrose availability, thereby linking growth and development to carbon status. The amounts of Tre6P in plant tissues are extremely low (between 10 and 1000 pmol-g⁻¹ FW in *Arabidopsis thaliana*), and plants contain several isomeric disaccharide monophosphates. HPLC coupled to tandem mass spectrometry (LC-MS/MS) is the most suitable technique for measuring Tre6P, offering both high sensitivity and high specificity.

Methods

Frozen plant material is extracted with chloroform/methanol(3:7, v/v) as in [1]. Each sample is spiked with a mix of [²H₂]Tre6P and other stable-isotope labelled internal standards (SIL-IS). Compounds are initially separated on an IonPac AS11-HC-4 with an alkaline eluent containing 5% methanol (ICS 5000+; Thermo Scientific Dionex). Detection and quantification is then achieved by tandem mass spectrometry (QTrap 6500; Sciex) [1,2].

Results and Conclusions

We achieve high specificity by simultaneous measurement of multiple product ions in the third quadrupole. The use of SIL-IS allows correction for ion suppression and other matrix effects in plant extracts which typically contain up to 25,000 different compounds [3]. In addition to Tre6P, we concurrently measure other phosphorylated intermediates, nucleotide sugars and organic acids, providing a comprehensive overview of central carbon metabolism.

Novel Aspect

This technique allows accurate determination of low abundance signal metabolites and metabolic intermediates in plants that cannot be measured using other methods.

References

1. Lunn JE, Feil R, Hendriks JH, Gibon Y, Morcuende R, Osuna D, Scheible WR, Carillo P, Hajirezaei MR, Stitt M. *Biochemical Journal* 397, 139–148, 2006.
2. Figueroa CM, Feil R, Ishihara H, Watanabe M, Kölling K, Krause U, Höhne M, Encke B, Plaxton WC, Zeeman SC, Li Z, Schulze WX, Hoefgen R, Stitt M, Lunn JE. *The Plant Journal* 85, 410–423, 2016.
3. Richard N Trethewey. *Encyclopedia of Life Sciences*, 1–8, 2007.

TP-186 / IDENTIFICATION OF NOVEL LONG CHAIN *N*-ACYLHOMOSERINE LACTONES OF CHAIN LENGTH C₂₀ FROM THE MARINE PHOTOTROPHIC BACTERIUM *RHODOVULUM SULFIDOPHILUM*

Naoko Arashida⁽¹⁾ - **Kazutaka Shimbo**⁽¹⁾ - **Takeshi Terada**⁽²⁾ - **Takuya Okimi**⁽²⁾ - **Yo Kikuchi**⁽³⁾ - **Shuhei Hashiro**⁽¹⁾ - **So Umekage**⁽²⁾ - **Hisashi Yasueda**⁽¹⁾

Ajinomoto Co., Inc., Institute for Innovation, Kawasaki (1) - Toyohashi University of Technology, Department of Environmental and Life Sciences, Toyohashi (2) - Waseda University, Graduate School of Science and Engineering, Tokyo (3)

Keywords: *N*-acylhomoserine lactone (AHL), quorum sensing (QS), C₂₀-HSL, precursor ion scan, high resolution mass spectrometry

Introduction

Gram negative bacterial quorum sensing (QS) is mainly regulated by an extracellularly produced *N*-acylhomoserine lactone (AHL). AHL consists of a lactone ring and an acyl chain, which generally varies from C₄ to C₁₈ in length. The marine photosynthetic bacterium *Rhodovulum sulfidophilum* produces

extracellular nucleic acids involved in its flocculation [1], and may be regulated through AHL signaling.

Methods

AHL extracts from cultured by *R. sulfidophilum* were obtained by TLC and analysed with liquid chromatography tandem mass spectrometry (LC-MS/MS) systems. In order to profile and detect AHLs, precursor ion scanning was performed using quadrupole-ion trap mass spectrometer. For both identification of the compounds and differential profiling analysis, high resolution accurate mass data was obtained by quadrupole-Orbitrap hybrid mass spectrometer.

Results

AHL fraction containing long acyl-chain homoserinlactone (HSL) exceeding C10 was obtained from the reverse-phase thin layer chromatography-fractionated cultured supernatant of the marine photosynthetic bacterium *R. sulfidophilum*. We performed precursor ion scanning analysis and detected two kinds of long chain AHLs with chain length C20. By fragmentation search analysis to detect compounds with a homoserine lactone ring moiety for data dependent acquisition, a minor AHL, presumed to be 3-OH-C18-homoserine lactone, was also found. Among the detected C20-HSLs, *N*-(3-Hydroxyicosanoyl)-L-homoserine lactone (3-OH-C20-HSL) was structurally identified and 3-OH-C20:1-HSL was strongly suggested (paper in submission).

Conclusions

We found that *R. sulfidophilum* produced two novel C20 long chain AHLs.

Novel Aspect

This is the first report on structural analysis of novel AHLs with the longest yet-described acyl chain length of C20.

References

1. Suzuki H, Daimon M, Awano T, Umekage S. et al. Appl Microbiol Biotechnol, 84, 349 (2009)

TP-187 / FATTY ACID β -OXIDATION FUELS TRICARBOXYLIC ACID CYCLE TO SUSTAIN CELL PROLIFERATION

Matteo Audano - Silvia Pedretti - Donatella Caruso - Nico Mitro - Emma De Fabiani

Dipartimento di Scienze Farmacologiche e Biomolecolari (DiSFeB), Università degli Studi di Milano, Milano, Italy

Keywords: Mitochondria, mtDNA, fluxomics, β -oxidation.

Introduction

Mitochondria play a central role in energy metabolism and possess their own genome (mtDNA). Mitochondria control cell cycle and differentiation where they orchestrate and finely tune all metabolic transitions. In line with this, mtDNA is vital for cell life since it encodes for 13 OXPHOS subunits, mt-tRNAs and mt-rRNAs. The goal of this project is to identify new metabolic patterns by which cells survive to an important metabolic deficit, such as lack of mtDNA.

Methods

For this study we took advantage of p0 cells (mtDNA partially depleted cells). By LC-MS/MS, we performed steady state metabolomics and metabolic tracing. Specifically, we labelled the cells with U-¹³C₆-glucose, U-¹³C₅-glutamine and U-¹³C₁₆-palmitate.

Results

Our results indicate that C2C12, Neuro2a, Hepa 1-6 and C3H/10T1/2 undifferentiated p0 cells have

different metabolic phenotypes compared to wild type control cells with normal mtDNA levels. Specifically, steady state metabolomics and metabolic tracing analysis indicate that mtDNA downregulation differentially impairs metabolic profile of p0 cells in a cell specific fashion. Specifically, metabolomic and fluxomic analyses indicate that β -oxidation fuels the TCA cycle to replenish the reducing cofactor NADPH through the activity of the malic enzyme 1 (ME1).

Conclusions

Our data indicate that mtDNA depletion to a similar extent, differentially affects the metabolic profile, depending on the cell type.

Novel Aspect

Our findings provide new insights into undisclosed metabolic pathways activated to face important metabolic deficits. Specifically, we directly link for the first time the β -oxidation of fatty acid to NADPH, aspartate and asparagine levels and cell proliferation.

TP-188 / DIAGNOSTIC PERFORMANCE OF METHYLMALONIC ACID TESTING BY LCMS/MS COMPARED TO HOLOTRANSCOBALAMIN BY IMMUNOASSAY FOR VITAMIN B₁₂ DEFICIENCY

Flair Bahadory - Phoebe Stanford - Rita A. Horvath

NSW Health Pathology, Prince of Wales Hospital, Department of Clinical Chemistry, Sydney, Australia

Keywords: Methylmalonic acid, LC-MS/MS, Holotranscobalamin, Immunoassay, Vitamin-B12 deficiency

Introduction

Measurement of total-B12 has been suggested to lack sensitivity to detect vitamin B₁₂ deficiency. Methylmalonic acid (MMA) and Holotranscobalamin (Holo-TC) are recommended to be more sensitive markers for B₁₂ deficiency. Our preliminary aim was to develop an LC-MS/MS method to measure MMA. We attempted to evaluate the diagnostic accuracy of HoloTC for the detection of B₁₂ deficiency based on the surrogate reference criterion of a raised MMA (>500nmol/L).

Methods

We employed a simple protein precipitation technique using acid, and a fast-chromatographic method using Ascentis RP-amide column and gradient mobile phase (3% - 90% v/v methanol with 0.1% formic acid) on Sciex 6500 QTRAP in negative mode, MRM of 117 - 73. Serum samples were tested for MMA (LC-MS/MS) and Holo-TC (immunoassay, Abbott-Architect / Roch-Cobas). Using ROC curve analysis and MMA as surrogate, the diagnostic accuracy of Holo-TC was evaluated.

Results

We successfully developed an LC-MS/MS method for MMA measurement with an analytical range of 0.14 - 86.4 μ mol/L. Serum samples, n=50 including 10 normal (>250pmol/L), 10 abnormal (<140pmol/L), 30 equivocal (140-250pmol/L) total B12 were tested for MMA and Holo-TC. Overall, we observed an increase in MMA (>500nmol/L) in 10% of patients, excluding patients with renal failure or known alternative reason for elevated MMA (1). Of this group, Holo-TC was low (<25 pmol/L) in 20% and equivocal (25-50pmol/L) in 80%. When MMA cut-off is 500nmol/L, the positive predictive value of Holo-TC<25pmol/L was 20%, and the negative predictive value of Holo-TC>50pmol/L was 100%. Compared to Abbott, results for Holo-TC by Roche was discordant in 4% of patients. All patients with low total B12 (<140pmol/L) had equivocal Holo-TC and normal MMA (<500nmol/L).

Conclusions

Both MMA and Holo-TC have limited abilities to predict B₁₂ deficiency. The use of these markers in combination, together with clinical details, may provide the most reliable representation of B₁₂ status.

Novel Aspect

To our knowledge, this is the first study that comparing several methods of measurement for detection of vitamin B₁₂ deficiency.

References

1. Loikas S, et al., Clin Chem Lab Med, 45, 197-201 (2007).

TP-189 / PYCMETRIX A FRAMEWORK FOR GAS CHROMATOGRAPHY - MASS SPECTROMETRY FOR CLINICAL DATA PRE-PROCESSING, ANALYSIS AND FEATURE EXTRACTION

Andrei Barcaru⁽¹⁾ - Rebecca Heiner-Fokkema⁽²⁾ - Péter Horvatovich⁽³⁾

University Medical Center of Groningen, Laboratory of medicine, Groningen (1) - University of Groningen, Center for Metabolic and Liver Diseases, Laboratory of medicine, Groningen (2) - University of Groningen, Interfaculty of Mass Spectrometry, Analytical Biochemistry group, Groningen (3)

Keywords: chemometrics, GC-MS data analysis, framework, python, Bayesian

Introduction

Gas chromatography coupled to mass spectrometry (GC-MS) is widely used in clinical data analysis due to several major advantages: speed, resolution and available databases. Several data preprocessing tools for GC-MS are available in R [1] and C [2]. However, free open source pipelines are less common in this field. Clinical data requires often times certain specifications and code editing thus emphasizing the need of open sourcing.

Methods

Pycmetrix is a framework build in python for preprocessing of GC-MS data and chemometrix data analysis. The package includes the possibility of deconvolution, quantification and interaction with a database via mysql.

Results

The framework was applied to GC-MS data processing of several patients diagnosed with MCADD. The quantification feature was validated using the calibration curves and showed very good results. A feature selection method was further developed for a fast and automated diagnosis. The data from 20 healthy patients and 20 MCADD were used to create a model in a lower feature space. Using the tool described herein, we were capable to accurately quantify and classify the data.

Conclusions

The pipeline proved to be efficient in working with multiple data sets in one run. The hardware requirements for the data analysis were reasonable, making the preprocessing and data analysis possible even on a low cost PC or laptop.

Novel Aspect

Pycmetrix includes the features not available in other python packages such as: effect size based feature selection, supervised projection pursuit, Bayesian hypothesis testing and other useful options.

References

1. Domingo-Almenara X., Brezmes J., Vinaixa M., Samino S., Ramirez N., Ramon-Krauel M., Lerin C., Díaz M., Ibáñez L., Correig X., Perera-Lluna A., Yanes O., Anal. Chem. 2016, 88, 9821-9829
2. Lommen A., Anal Chem. 2009, 81, 3079-86

TP-190 / THE UTILIZATION OF TANDEM MASS SPECTROMETRY TECHNIQUES IN DIAGNOSIS OF PATIENTS WITH HOMOCYSTINURIAS.

Josef Bártil - Jakub Krijt - Petr Chrastina - Jakub Hodík - Markéta Paulová - Renata Pinkasová - Jitka Sokolová - Viktor Kožich - Karolína Peškova

General University Hospital in Prague and 1st Faculty of Medicine, Charles University in Prague, Clinic of Pediatrics and Adolescent Medicine, Prague

Keywords: tandem mass spectrometry, homocystinuria

Introduction

Homocystinurias (HCYU) are rare genetic disorders characterized by accumulation of homocysteine (Hcy) in the body due to defects in a) the remethylation (RMD) pathway or b) trans sulfuration pathway (cystathionine β -synthase deficiency -CBSD). Metabolism of Hcy depends also on supply, transport and intracellular metabolism of B vitamins. We developed tandem mass spectrometry methods enabling differential diagnosis of various types of HCYU.

Methods

FIA-MS/MS newborn screening (NBS) method for amino acids (AA) in dried blood spot (DBS) and LC-MS/MS determinations of: 1) Hcy, methionine (Met) and methylmalonate (MMA) in DBS. 2) Cystathionine and other AA in plasma after derivatization with alkyl chloroformates. 3) S-adenosyl methionine (SAM) and S-adenosyl homocysteine (SAH) in plasma. 4) CBS activity in plasma using stable isotope substrate and determination of the deuterated product cystathionine D2.

Results

The established reference ranges were as follows: Hcy (DBS): < 12 $\mu\text{mol/L}$; Met (DBS): 7-33 $\mu\text{mol/L}$; MMA (DBS): < 5 $\mu\text{mol/L}$; Hcy (plasma): 3.5-15 $\mu\text{mol/L}$; cystathionine in plasma (P): 80-1000 $\mu\text{mol/L}$; Met (P): 12-40 $\mu\text{mol/L}$; SAM (P): 50-160 nmol/L; SAH (P): 26-80 nmol/L and CBS activity: 100-1000 nmol/hour/L plasma. By introducing new LC-MS/MS methods we diagnosed patients with CBSD presented by accumulation of Hcy and Met and decreased metabolites below the CBS block (*i.e.* cystathionine and cysteine), RMD including methylene tetrahydrofolate reductase deficiency and intracellular methylcobalamin synthesis defects characterized by decreased production of Met and SAM, accumulation of Hcy and an increased flux of metabolites through the trans sulfuration pathway and B vitamins deficiencies presenting with elevated Hcy and MMA. Within NBS program between July 2016 and December 2017, we analyzed samples from 128531 newborns and detected 3 patients with maternal vitamin B12 deficiency and 1 patient with CBS deficiency.

Conclusions

We developed LC-MS/MS methods approach for diagnosing patients with HCYU enabling us to discriminate between primary enzyme deficiencies and vitamin B deficiencies and distinguish between defects in the RMD and trans sulfuration pathway of Hcy. methods employing DBS are also used in NBS program and have the clinical utility for monitoring of therapy in patients with hyperhomocysteinemias due to ease of collecting and shipping DBS samples by the patients.

Novel Aspect

The LC-MS/MS methods help to enhance efficacy in diagnosis procedure for hyperhomocysteinemias and simplify monitoring of patients under treatment.

Acknowledgement

This work was supported by the grant Nr. 16-30384A from the Czech Health Research Council and RVO VFN64165 from the Ministry of Health of the Czech Republic.

TP-191 / CESI-MS - A SENSITIVE AND VERSATILE APPROACH FOR METABOLOMICS

Bryan Fonslow⁽¹⁾ - ***Esme Candish***⁽¹⁾ - ***Jim Thorn***⁽²⁾ - ***Walter Belloni***⁽³⁾ - ***Stephen Lock***⁽²⁾

SCIEX, Framingham, MA (1) - SCIEX, Warrington, United Kingdom (2)

Keywords: CE-MS, anionic, cationic, metabolites, CESI

Introduction

The analysis of the metabolites can present an analytical challenge, particularly those involved in the central carbon metabolism and neurotransmission. Often, structural isomers can't be differentiated by MS or MS/MS alone, making a separation crucial. Unfortunately, the LC separations prove difficult for some metabolites; some hydrophilic analytes are poorly retained by RP-LC, while column-to-column reproducibility can sometimes be more challenging with HILIC separations.

Methods

In contrast, capillary electrophoresis (CE) is well suited towards metabolomics analyses. CESI-MS integrates CE and electrospray ionization (ESI) into a single device and this combines the benefits of a high-resolution separation with the increased MS sensitivities that result from the ultra-low flow rates. In this case, we demonstrate CESI-MS for the analysis of these challenging metabolites.

Results

A single protocol was established for the separation and detection of both cationic and anionic metabolites. Electrophoretic separations were performed using 30 kV to generate a field strength of 333 V/cm and normal and reversed CE polarities were employed for cationic and anionic metabolites respectively. MS compatible background electrolytes (BGEs) facilitated highly efficient separations, achieved in less than 30 mins. Enhanced sensitivity was achieved using transient isotachopheresis with CESI directly coupled to either a SCIEX TripleTOF® 6600 or QTRAP® 6500+ system. Cationic metabolites were characterized using normal polarity with positive ESI while the anionic metabolites were analyzed from the same sample by simply switching the polarity of the CE and MS.

Conclusions

The CESI-MS technology proved powerful for the targeted analysis of the central carbon metabolism, particularly the small organic acids of the tricarboxylic acid (TCA) cycle and the isobaric phosphorylated sugars of the pentose and glycolysis pathways. In addition, CESI-MS was used to detect metabolites involved in neurotransmission such as those in the tryptophan pathway. Despite an injection volume of less than 50 nL, low nanoMolar concentrations were easily detected.

TP-192 / ASSESSMENT OF KIDNEY GRAFT QUALITY BY SPME-(LC)-HRMS: GOING BEYOND CLINICAL ROUTINE

Iga Stryjak⁽¹⁾ - ***Matyas Hamar***⁽²⁾ - ***Markus Selzner***⁽²⁾ - ***Barbara Bojko***⁽¹⁾

Nicolaus Copernicus University in Torun, Department of Pharmacodynamics and Molecular Pharmacology, Faculty of Pharmacy Collegium Medicum in Bydgoszcz, Bydgoszcz (1) - Toronto General Hospital, University Health Network, Multi Organ Transplant Program, Department of Surgery, Toronto (2)

Keywords: solid phase microextraction, high resolution mass spectrometry, metabolomics, transplantation

Introduction

Since successful introduction of immunosuppressants, organ transplantation is a routine procedure, which saves lives of patients with end-stage organ failure worldwide. The demand for grafts is much higher than

the available pool of organs, therefore various strategies are explored to address this issue. The one presented here is metabolic profiling of organ with low-invasive method based on solid phase microextraction and high resolution mass spectrometry.

Methods

The study of kidney grafts were performed on pig model. In one group organs were subjected to warm ischemia and in the second group kidneys were non-ischemic. Solid phase microextraction probes of ca. 0.2 mm diameter and 7 mm of coating length were selected for untargeted extraction of metabolites. Sampling was performed before organ harvesting, during preservation with normothermic perfusion and after reperfusion. Analysis were done on UPLC-Q-Exactive Focus.

Results

The proposed analytical strategy enabled metabolomics of living organ with very low-invasiveness and tissue biopsy-free sampling. The results of LC-MS analysis revealed significant differences in metabolome of kidney obtained from donors after cardiac death (DCD, ischemic organs) and heart beating donors (no ischemia). Comparison of the data from organs subjected to warm ischemia showed increased level of lysyl residue oxidation products. The effect was most pronounced after reperfusion indicating increased risk of reperfusion injury in DCD. On the other hand, monitoring of changes during normothermic perfusion showed that level of some metabolites inducing oxidative stress peaking during ischemia drops to initial values suggesting that the used preservation strategy can successfully reverse some adverse reactions initiated by ischemia. This finding is supported by putative identification of the nephroprotective compounds reported earlier in the literature as the ones possessing anti-ischemia-reperfusion injury activity.

Conclusions

The results obtained indicate that novel analytical strategy proposed for assessment of graft quality carries a great potential for future improvement of pre-transplant diagnostics. It offers a possibility of selection of panels of biomarkers, which after required validation, could be used for fast on-site screening with SPME-MS platform in the future.

Novel Aspect

The applicability of SPME-HRMS-based strategy for graft assessment was proven using a porcine model of kidney transplant.

TP-193 / MULTILAYER OMICS STRATEGY FOR TRACING THE METABOLIC ROUTES EXPLOITED BY KRAS ONCOGENETIC MUTATION IN NSCLC

Laura Brunelli⁽¹⁾ - ***Elisa Caiola***⁽²⁾ - ***Francesca Falcetta***⁽²⁾ - ***Silvia Giordano***⁽¹⁾ - ***Mirko Marabese***⁽²⁾ - ***Massimo Brogini***⁽²⁾ - ***Roberta Pastorelli***⁽¹⁾

IRCCS-Istituto di Ricerche Farmacologiche Mario Negri, Environmental Health Sciences, Milano (1) - IRCCS-Istituto di Ricerche Farmacologiche Mario Negri, Oncology, Milano (2)

Keywords: NSCLC, metabolism, co-occurring genetic lesions, mass spectrometry

Introduction

NSCLC is a heterogeneous disease, with multiple oncogenic mutations (*e.g.* KRAS, LKB1). NSCLC patients have a variety of clinical courses, which are partly due to differences in co-occurring molecular events [1]. Since both KRAS mutation and LKB1 loss have an impact on cellular metabolism, it is fundamental to discern the metabolic effects induced by single genetic events from those induced by their co-occurrence, for novel therapeutic interventions

Methods

We knocked-out the LKB1 gene in well-characterized NSCLC cell clones harbouring KRAS WT or mutant G12C proteins. We obtained an isogenic system in which KRAS mutation and LKB1 loss were individually or concomitantly present. The effects of these genetic makeups on cell metabolism were investigated in these NSCLC cells by means of an integrated survey of proteomics, metabolomics, ¹³C-Glucose/Glutamine flux analysis coupled with functional in-vitro strategies.

Results

Using untargeted (FIA-LTQ Orbitrap) and targeted (p180KIT, HPLC-Triple Quad 5500) metabolomics approaches we demonstrated how isogenic NSCLC cells harbouring both KRASG12C/LKB1loss had an impact on the same metabolic pathways (glycolysis, TCA, urea cycle) used by their single lesion counterparts, but were able to exploit these metabolic routes through a heightened metabolite production. The enhanced metabolite production triggered by the co-occurrence was not reflected in a different enzymatic asset in terms of both numbers and magnitude of enzymatic changes. Furthermore, metabolic flux analysis showed how the metabolic dynamics of glycolysis and TCA cycle differed between clones harbouring the dual lesion from those with the single lesions. The enhanced metabolic activity of cells with both genetic lesions rendered the viability of these cells susceptible to energetic stress caused by both glycolysis inhibition and by nutrients limitation (glucose/glutamine) than those harboring single oncogenic lesions.

Conclusions

We demonstrate how the metabolic rewiring observed in presence of KRAS/LKB1 co-occurring lesions, whilst failing to improve cellular fitness, was accompanied by reduced cell growth rates compared to those of wild type or single lesion cells. Furthermore, the cells with both lesions demonstrated greater susceptibility towards survival impedance caused by 2-DG treatment or glucose/glutamine limitation.

Novel Aspect

This observation raises the prospect that energy stress may affect NSCLC cells harbouring co-occurring lesions, which may render them more susceptible to cytotoxic drugs

References

1. Bhattacharya S. Socinski M. A., Burnset T. ClinTransl Med. 4, 35; 2015

TP-194 / CONVERSION OF COENZYME A THIOESTERS BY NATIVE CHEMICAL LIGATION – A NEW APPROACH FOR CHARACTERIZING AN IMPORTANT CLASS OF BIOMOLECULES

Nevenka Cakić - Heinz Wilkes

Institute for Chemistry and Biology of the Marine Environment (ICBM), University of Oldenburg, Oldenburg

Keywords: Coenzyme A thioesters, native chemical ligation, fragmentation patterns, chiral separation

Introduction

Acyl-coenzyme A esters (acyl-CoAs) are very important biomolecules which play a central role in different metabolic processes, including β -oxidation of fatty acids, amino acid degradation and the citric acid cycle. Although numerous techniques for the analysis of acyl-CoAs have been developed, a straightforward approach for the unambiguous characterization of isomeric acyl-CoAs is still lacking.

Methods

Transformation of CoA-thioesters has been achieved using the native chemical ligation (NCL) reaction [1]. Different acyl-CoAs were incubated with L-cysteine and sodium 2-mercaptoethanesulfonate (MESNa) in 0.05M phosphate buffer (pH 7.4) at 70 °C. Products were analyzed without further modification and after the cleavage of the disulphide bond by a reductive agent with LC-ESI-TOF in positive mode.

Results

In this research the replacement of coenzyme A from CoA-thioesters with L-cysteine by NCL is demonstrated. Initially, transthioesterification between acyl-CoAs and L-cysteine occurs to form products that further undergo an irreversible amide formation. Under the reaction conditions applied the thiol-containing amide is oxidized to a disulphide with MESNa, which results in higher water solubility. All products detected showed mass spectrometric fragmentation patterns corresponding to cleavage of the amide and disulphide bonds, which allow their straightforward identification in mixtures. After reduction of the disulphides all products showed a fragmentation pattern characteristic of cysteine derivatives. Additionally, we performed the NCL reaction on the racemate of 2-methylhexanoyl-CoA and successfully separated diastereomeric products on an achiral column.

Conclusions

The methodology described enables detailed structural characterisation of the acyl-residues in acyl-CoAs by generating transformation products through the NCL reaction which deliver more informative mass spectra. Furthermore, it provides a possibility for separation and assignment of configuration of alpha-substituted chiral compounds on an achiral column using L-cysteine as a chiral auxiliary.

Novel Aspect

Selective conversion of CoA thioesters to the cysteine containing products is a novel aspect for characterization of an important class of biomolecules.

References

1. Kulkarni S. S., Sayers J., Premdjee B. and Payne R. J., Nature Reviews Chemistry 2, 0122 (2018).

TP-195 / ENHANCED METABOLITE IDENTIFICATION USING AN ORBITRAP TRIBRID MASS SPECTROMETER

Kate Comstock⁽¹⁾ - *Shuguang Ma*⁽²⁾ - *Seema Sharma*⁽¹⁾ - *Helene Cardasis*⁽¹⁾ - *Yan Chen*⁽¹⁾

Thermo Fisher Scientific, San Jose (1) - Genentech Inc., South San Francisco (2)

Introduction

High resolution mass spectrometry (HRMS) is an essential tool for metabolite identification. HRMS with sophisticated data acquisition features can provide the critical information for metabolite structure elucidation and alleviate the difficulties of matrix complexity.

Here we present a study for metabolite ID using Orbitrap Tribrid MS with data acquisition features designed to significantly improve small molecule structure analysis, which include automatic background ions subtraction, inclusion/exclusion lists generation, and real-time decision-making to optimize MSn data acquisition quality and speed.

The high quality information-rich MSn data generated using this improved HRMS was processed using structure analysis and data mining software "Mass Frontier" and "Compound Discoverer".

Methods

Amprenavir, Bosentan, Lopinavir, and Ritonavir were selected as model compounds. The compounds, at concentration of (5 µM), were incubated with human liver microsomes (HLM) in the presence of GSH and UDPGA. Incubations without drug or without NADPH were chosen purposely to test the background subtraction feature.

The LCMS analyses were performed on Vanquish Flex Binary UPLC system with DAD detector coupled to an Orbitrap Tribrid Mass Spectrometer. Mobile phases were composed of: A: H₂O/0.1% formic acid, and B: ACN/0.1% formic acid with gradient on Hypersil column (2.1×100 mm 1.9 µm).

High resolution full scan and MSn data were collected in a data-dependent fashion with polarity switching. The data acquisition method was incorporated with features optimized for small molecule analysis.

Preliminary Data

These HIV drugs were selected because of their extensive metabolism. Using the Orbitrap Tribrid Mass Spectrometer's special features for small molecule data acquisition: the automatic background subtraction, inclusion/exclusion list generation, and real-time decision-making optimization of MS_n spectra acquisition, the MS provides a wealth of information for metabolite identification. The automatic background subtraction feature is especially useful for identification of metabolites present in a complicated biologic matrix.

The new and enhanced features that use the high quality MS/MS fragments to trigger the MS_n acquisition significantly improved MS_n data quality, because multi-stage MS/MS (MS_n) information is essential for confident metabolite structure elucidation. For example, HRAM full scan data show that Ritonavir metabolites include four mono-oxidation products, two *N*-dealkylation products, and four di-oxidation products. In order to determine the structure and the sites of biotransformation of these metabolites, high quality MS/MS and higher order MS/MS spectra were obtained using the novel data collection features of Orbitrap Tribrid MS.

Major and minor metabolites, both phase I and II, of these model compounds formed in the HLM microsomal incubations were quickly detected by processing the HRAM and higher level MS/MS (MS_n) spectra of metabolites using Mass Frontier and Compound Discoverer.

This workflow offers significant improvements in speed and confidence of routine drug metabolite identification and other small molecule structure characterization applications.

Novel Aspect

Small molecule specific acquisition features of Orbitrap Tribrid MS improve MS_n data quality and increase confidence for metabolite identification.

TP-196 / TARGETED APPROACH FOR DETERMINATION OF METABOLITES OF FLAVONOIDS AND ISOFLAVONOIDS IN BIOLOGICAL SAMPLES

Lucia Chrenková⁽¹⁾ - *Veronika Pilařová*⁽¹⁾ - *Lenka Applová*⁽²⁾ - *Iveta Najmanová*⁽²⁾ - *Přemysl Mladěnka*⁽³⁾ - *Lucie Nováková*⁽¹⁾

Charles University, Faculty of Pharmacy, Department of Analytical Chemistry, Hradec Králové (1) - Charles University, Faculty of Pharmacy, Department of Pharmacology and Toxicology, Hradec Králové (2) - Charles University, Faculty of Pharmacy, Department of Biological and Medical Sciences, Hradec Králové (3)

Keywords: quercetin, UHPLC-MS/MS, MEPS, μ -SPE-PT, isoflavonoids

Introduction

Rutin is metabolized into two groups of metabolites: flavonoid type with higher molecular weight (quercetin, quercetin-3-O-glucuronide, tamarixetin and isorhamnetin) and polar weak phenolic acids (4-methylcatechol, phloroglucinol, 3,4-dihydroxyphenylacetic acid, homovanilic acid, 3-hydroxyphenylacetic acid, 3-(3-hydroxyphenyl)propionic acid). Tectoridin and tectorigenin are isoflavonoids.

Methods

The UHPLC method with mass spectrometry (MS/MS) detection was optimized for the separation of all 13 analytes. For the pre-treatment of rat plasma samples, microscale extraction approaches including μ -SPE-PT (micro-solid phase extraction in pipette tips) and MEPS (microextraction by packed sorbent), but also non-selective PP (protein precipitation) were used.

Results

The metabolites with significantly different physicochemical properties make the simultaneous determination quite challenging.

Considering the critical pair of isomers (tamarixetine and isorhamnetin), BEH Shield RP C18 (2.1 x 100 mm; 1.7 μ m) column and gradient elution with methanol and 0.1% aqueous solution of formic acid offered the best selectivity. The Xevo TQ-XS mass spectrometry was employed and the conditions for all compounds were also thoroughly optimized. Due the instability of some phenolic acids (homovanilic acid, 3,4-dihydroxyphenylacetic acid) in ion source, the Xevo TQ-S mass spectrometer detector was replaced with Xevo TQ-XS. This mass spectrometer contains different type of ion optics more suitable for targeted analytes. All compounds were analyzed in the mode SRM (selective reaction monitoring), while SIM was often needed at other platforms.

Conclusions

The UHPLC-MS/MS method was optimized and validated in the terms of sensitivity, linearity and matrix effects. Using Xevo TQ-XS mass spectrometer, all analytes could be quantified with linearity range between 0.1 – 1000 ng/mL.

Novel Aspect

Due the instability of some phenolic acid analytes in ion optics, the Xevo TQ-S detector was replaced with Xevo TQ-XS.

This work was supported by the STARSS project (Reg. No. CZ.02.1.01/0.0/0.0/15_003/0000465) co-funded by ERDF, by SVV 260412 and by GACR 301/17/05409S.

TP-197 / EXPLORING SERUM TARGETED METABOLOMICS MASS SPECTROMETRY PROFILE ANALYSIS TO IMPLEMENT RISK CRITERIA FOR FIRST-DEGREE RELATIVES (FDR) OF GASTRIC CANCER (GC) PATIENTS

Giuseppe Corona ⁽¹⁾ - ***Renato Cannizzaro*** ⁽²⁾ - ***Gianmaria Miolo*** ⁽³⁾ - ***Agostino Steffan*** ⁽⁴⁾ - ***Valli De Re*** ⁽⁴⁾

IRCCS- National Cancer Institute, Immunopathology and Cancer Biomarkers, Aviano (PN) (1) - National Cancer Institute, Oncological Gastroenterology, Aviano (2) - National Cancer Institute, Medical and Preventive Oncology, Aviano (PN) (3) - National Cancer Institute, Immunopathology and Cancer Biomarkers, Aviano (PN) (4)

Keywords: Metabolomics; Gastric cancer; first degree relatives; Biomarkers; Early diagnosis; introduce

Introduction

A positive family history is a strong and consistently reported risk factor for gastric cancer(GC) [1]. The main purpose of this study was to investigate serum metabolomic profiles to find additional biomarkers that could be integrated with serum pepsinogens and gastin test to improve the diagnosis of GC and the selection of first-degree relatives (FDR) at higher risk of GC development.

Methods

The targeted serum metabolomic profile included 188 metabolites, covering amino acids, biogenic amines, acylcarnitines, phosphatidylcholines, sphingomyelins and hexoses. Metabolite profiles were performed by liquid chromatography and flow injection tandem mass spectrometry analysis. Differential metabolomic signatures among the FDR (n=54)and GC (n=71) groups were investigated by univariate and multivariate analysis.

Results

Forty metabolites mainly belong mainly the phospholipids and acylcarnitines classes were significantly altered between FDR and GC patients. Nine metabolites resulting from the training set were further confirmed in the validation set. Compared with FDR, GC patients were characterized by lower levels of hydroxylated sphingomyelins (SM(OH)22:1, SM(OH)22:2, SM(OH) 24:1) and phosphatidylcholines (PC ae 40:1, PC ae 42:2, PC ae 42:3) and by higher levels of acylcarnitines derivatives (C2, C16, C18:1). The

sensitivity and sensibility of these metabolomic biomarkers to distinguish potential cancer conditions was 73.47% and 83.78% respectively with AUC of ROC curve of 0.811 that improves to 0.90 when integrated with age and the serum level of serum pepsinogen II.

Conclusions

The results of the study for the first-time described serum metabolomic profiles that discriminate GC patients from FDR sharing the same environment and a similar genetic background. As compared with FDR, the GC patients showed specific serum metabolomic signatures characterized by an increase in specific acylcarnitines and a decrease in a distinctive subclass of sphingolipids.

Novel Aspect

The study underlines the role of the use of the individual's serum tandem mass spectrometry metabolomics profile to complement the triage of FDR at higher risk of GC development

References

1. Torre, L. A. ; Siegel, R. L. ; Ward, E. M. ; Jemal, A. Cancer Epidemiol. Biomarker Prev., 25, 16–27,2016.

TP-198 / PHYTOCHEMICAL SCREENING OF FAGIOLI BIANCHI DI ROTONDA BEANS PDO (PHASEOLUS VULGARIS L.) AND EVALUATION OF THEIR BIOLOGICAL ACTIVITY FOR FUNCTIONAL FOOD PROPERTIES

Donatella Coviello⁽¹⁾ - ***Alberto Onzo***⁽¹⁾ - ***Giuliana Bianco***⁽¹⁾ - ***Raffaella Pascale***⁽²⁾ - ***Maria Cristina Lafiosca***⁽¹⁾ - ***Philippe Schmitt-Kopplin***⁽³⁾ - ***Patrizia Iannece***⁽⁴⁾ - ***Carmine Gaeta***⁽⁴⁾

Dipartimento di Scienze, Università degli Studi della Basilicata, Potenza (1) - Scuola di Ingegneria, Università degli Studi della Basilicata, Potenza (2) - Analytical BioGeoChemistry, Helmholtz Zentrum Munchen, Neuherberg (3) - Dipartimento di Chimica, Università degli Studi di Salerno, Fisciano (4)

Keywords: Fagioli Bianchi di Rotonda, FT-ICR-MS, Phytochemical Screening, Hypoglycemic Activity, Metabolomics

Introduction

The aim of this study was to evaluate the metabolome of Fagioli Bianchi di Rotonda beans (*Phaseolus Vulgaris* L.), cultivated in Basilicata (Southern Italy), labelled with DOP mark quality, in order to promote their nutraceutical properties. In this study, a metabolome untarget analysis by HRMS was performed.

Methods

High-resolution mass spectra were acquired on a Bruker (BrukerDaltonik GmbH, Bremen, Germany) solarix XR Fourier transform ion cyclotron resonance mass spectrometer (FT-ICR-MS) equipped with a 7 Tesla superconducting magnet and an ESI source in the negative ionization mode. The nutritional compositions including starch, fat, protein and phytochemicals were analyzed.

Results

A shotgun approach, based on direct infusion negative-ion ESI ultrahigh resolution mass spectrometry (FT-ICR-MS), was employed for the rapid analysis of metabolites occurring in all extracts of Fagioli Bianchi di Rotonda beans. Such non-targeted analysis generates a tremendous amount of data and requires visualization strategies to convert lists of accurate m/z values into metabolomic context. Visualization strategy using van Krevelen diagram was adopted. This diagram displays the hydrogen/carbon (H/C) vs. oxygen/carbon (O/C) ratios of these elemental formulas and provide a qualitative description of the molecular complexity of Fagioli Bianchi di Rotonda data. This plot enabled the localization of chemical species according to class metabolites, as carbohydrates, glycosylated compounds, peptides, polyphenols and saponins. The antioxidant activity and α -glucosidase inhibition activity of extracts were evaluated. All the samples possessed strong ABTS free-radical-scavenging capacity and α -glucosidase inhibition activity.

Conclusions

In this study, an identification of different classes of secondary metabolites in Fagioli Bianchi di Rotonda, ecotypes Bianco and Tondino, was carried out. Some of these compounds are known to show bioactivities interesting for human health. results of biological activity assays highlighted the importance of the use of such food product for nutraceutical applications.

Novel Aspect

This study shed some light on metabolic profiles of Fagioli Bianchi di Rotonda and their nutraceutical properties. results can be used for comparative studies.

References

1. Pascale R., Bianco G., Cataldi T. R. I., Schmitt-Kopplin P., Bosco F., Vignola L., Uhl J., Lucio M., Milella L., Food Chemistry 242 (2018) 497–504.
2. Bianco G., Pascale R., Carbone C. F., Acquavia M. A., Cataldi T. R. I., Schmitt-Kopplin P., Buchicchio A., Russo D., Milella L., Analytical and Bioanalytical Chemistry 410 (2018), 1561–1569.

TP-199 / TWO-DIMENSIONAL GAS CHROMATOGRAPHY AS A POWERFUL TOOL FOR DISCRIMINATING VOCs IN MONOVARIETAL AND COMMERCIAL EXTRA VIRGIN OLIVE OIL

Alessio Da Ros⁽¹⁾ - **Silvia Carlin**⁽¹⁾ - **Domenico Masuero**⁽¹⁾ - **Samantha Riccadonna**⁽¹⁾ - **Fulvio Mattivi**^(1,2) - **Nadia Mulinacci**⁽³⁾ - **Urska Vrhovsek**⁽¹⁾

Fondazione Edmund Mach, Food quality and Nutrition, San Michele all'Adige (TN) (1) - University of Trento, San Michele all'Adige (TN) (2) - NEUROFARBA, University of Florence, Florence (3)

Keywords: GCxGC-TOF-MS, EVOO, VOCs, oils quality, LOX pathways, terpenes.

Introduction

The unique and special flavour of Extra Virgin Olive Oil (EVOO) can be attributed to VOCs such as the C₆ and C₅ aldehydes, which are responsible for green aromatic notes, but also terpenes and sesquiterpenes, which tend to generate greater aromatic complexity in quality oils. The identification of the compounds responsible for flavour or off-flavour attributes is considered one of the keys for EVOOs quality control and a good indicator of olive oil quality changes.

Methods

Sixteen monovarietal EVOOs produced from different cultivars (purchased directly from producers) were compared with 25 blended commercial EVOOs (purchased in national grocery stores).

EVOO VOCs were analysed with HS-SPME-GCxGC-TOF-MS [1]. A total of 1254 putative compounds were detected, with 123 being significantly different in monovarietal and commercial oils. Compound identity was confirmed with Kovats retention index and by comparing spectra compounds to NIST 2.0, Wiley 8 FFNSC 2 (Chromaleont, Messina, Italy).

Results

Principal component analysis (PCA) showed clear separation between the two classes of samples. Specifically, commercial samples were more influenced by compounds such as alcohols, acids (acetic acid and formic acid) and other acetates frequently responsible for certain unpleasant notes. The monovarietal samples were more influenced by compounds deriving from LOX pathways (aldehydes C₆ and C₅), as well as terpenes and sesquiterpenes, such as valencene, copaene and ocimene for example, which generally generate pleasant wood, green and lemon notes. Furthermore, differences in the composition of monovarietal samples were evaluated. Two varieties (Casaliva and Tonda Iblea) were clearly distinguished from the rest of the monovarietal samples, underlining that the cultivar and pedoclimatic conditions can

influence VOC composition.

The preliminary results show that HS-SPME-GCxGC-TOF-MS is a fast and very powerful tool for revealing differences between the aroma of different classes of EVOO products.

Conclusions

The present investigation demonstrates that volatile compounds can be used as reliable indicators of Extra Virgin Olive Oil quality. In detail, between commercial and monovarietal oils, significant qualitative differences were found in a large number of volatile compounds, highlighting the possibility to use terpenes, sesquiterpenes and aldehydes how markers for quality olive oil aroma.

Novel Aspect

Two-dimensional gas chromatography is a novel method to discriminate EVOOs quality and aroma. Its applicability can reveal the differences of a large number of volatile compounds, analyzing the macro and micro differences in this matrix.

References

1. G. Purcaro, C. Cordero, E. Liberto, C. Bicchi, L. S. Conte, Toward a definition of blueprint of virgin olive oil by comprehensive two-dimensional gas chromatography, *J. Chromatogr. A* 1334, 101-111 (2014).

Acknowledgements

The research was supported by the Autonomous Province of Trento (Italy) and we thank Agririva – Agraria Riva del Garda for kindly providing the samples.

TP-200 / A TARGETED METABOLOMICS ASSAY TO INVESTIGATE STEROID HORMONES IN WOMEN DURING PREGNANCY AND EARLY POSTPARTUM PERIOD

Serena Di Palma - Endre Laczko

Functional Genomics Center Zurich, University of Zurich, ETH, Zurich

Keywords: Steroid hormones, saliva, pregnancy, targeted metabolomics, SRM

Introduction

In this study, we want to set up a targeted assay to measure the concentration of six steroids from salivary samples of pregnant women, beginning from the first trimester of pregnancy to a few weeks postpartum.

Saliva has potential clinical applications in maternal monitoring and represents a less invasive alternative to serum or plasma for the measurement of steroid hormones. However, little information is available regarding the range of salivary steroid concentrations during pregnancy [1, 2].

Steroid hormones in saliva derive primarily from free steroid present in the general circulation and their concentrations can be regarded as the bioavailable fraction. Our aim is to establish an accurate MS-based assay for quantification of six steroids in saliva, reaching great sensitivity and precision.

Methods

Steroid hormones are extracted from the biological matrix employing a methanol extraction protocol. After reconstitution in an aqueous buffer containing an internal standard, they are analyzed using a targeted approach by SRM technology on a nanoLC-MS system (Thermo TSQ-Quantiva coupled to a Waters M-class UPLC). Our targets include cortisol, cortisone, estradiol-17 β , estriol, progesterone, dehydroepiandrosterone (DHEA), and testosterone. A minimum of three transitions are selected per each target compound, with optimized collision energy, to ensure high selectivity and great sensitivity. The LC is set up in RP mode, with a C18 column (200 μ m x 5cm), a flow rate of 2 μ L/min and a gradient of 10 minutes. Calibration curves are included in the analysis for accurate quantification. The data will be investigated using a combination of Xcalibur Quan Browser and TraceFinder.

Results

We have fully developed and optimized an SRM assay to include our six steroid targets in the analysis. For the purpose, commercially available reference compounds were directly infused into the MS to select the best precursor ion and three to four fragment ions, optimizing collision energy (among other parameters) per each fragment. Also, the LC gradient was optimized to have each peak, corresponding to each compound, chromatographically well resolved.

Calibration curves containing all our compounds simultaneously, plus an internal standard to normalize for instrument variations, were built ranging from 0.4 to 400 nM. The curves showed good linearity and, for most of the compounds, the LLOQ was below 1 nM.

A test phase was conducted on saliva samples (100 µL) collected from pregnant women upon awakening at different stage of pregnancy. These test samples were employed to refine the extraction procedure and evaluate sample recovery.

We will expand the analysis to hundreds of saliva samples collected from pregnant women participants, monitoring the full range of steroid hormone concentration across pregnancy. After completing the acquisition, we will build quantitative plots of our targeted steroid compounds per each woman.

Conclusions

Targeted metabolomics employing SRM assays allows us to investigate a broad panel of steroid hormones in biological fluids, such as saliva sample. The study will shed light on the course of the steroid concentrations in saliva across pregnancy and into the early postpartum period.

Novel Aspect

Our MS assay will represent a fast and reliable tool to measure steroid concentrations in the saliva of pregnant women.

References

1. Nierop, A., Bratsikas, A., Klinkenberg, A., Nater, U. M., Zimmermann, R., Ehlert, U. *Journal of Clinical Endocrinology and Metabolism*, 91, 1329–1335 (2006).
2. Harris, A., Seckl, J. *Hormones and Behavior*, 59, 279–289 (2011).

TP-201 / ECOMETABOLOMICS OF ROOT EXUDATES IN GERMAN GRASSLAND COMMUNITIES

Sophie Dietz - Stefanie Döll - Karin Gorzolka - Dierk Scheel

Leibniz Institute of Plant Biochemistry, Stress and Developmental Biology, Halle-Saale

Keywords: Life Science, Metabolomics, LC-MS, untargeted profiling

Introduction

Rhizosphere: the belowground black-box.

The plant root is embedded in a multiple organism network: the rhizosphere. For reactions towards different abiotic stresses [1, 2] or to interact with the microbial soil communities and neighboring plants [3-7], roots release metabolites called exudates. Until now, there is little known about these chemicals in natural habitats. Therefore, we applied non-targeted metabolomics combined with ecology to analyze exudate patterns in natural grassland communities.

Methods

Exudate analysis under field condition.

Ten common European grassland species were grown under field conditions, exudates were directly collected from roots and subjected to non-targeted metabolite profiling by UPLC coupled to ESI-Q-TOF mass spectrometry and the combination of different software. Selected compounds were chemically classified on the level of metabolite families by the use of MetFamily. The method for exudate isolation was

proven to be appropriate by control experiments.

Results

Exudates: species specific and environmentally influenced.

Using our approach, more than 2800 compounds from samples of field grown plants were annotated. Different statistical tools revealed the release of semi-polar metabolites in a growth form dependent manner with species-specific compound patterns. This is more distinct in forbs than in grasses, which showed a less differentiated metabolite composition. In order to identify the main factors responsible for the exudate patterns, the species itself, plant functional traits, neighboring plants and environmental conditions were considered. In case of the five forbs, target species identity was determined as the main driver for exudate composition, while environmental conditions are decisive in case of the five grasses. Nevertheless, there are no specific metabolite classes released by single species. Rather, most of the detected specific metabolites could be putatively classified as members of the typically defense related metabolite families of flavonoids, coumarins, isoprenoids and phenylpropanoids, but also sugar related metabolites.

Conclusions

MS exudate analysis under field condition is feasible.

Our analysis demonstrates the possibility to obtain and analyze root exudates in field experiments. The species specificity of root exudate pattern appears either to be driven by interspecific differences or by the environmental conditions, depending on the investigated growth forms. The putative classification on the level of metabolite families delivered a broad spectrum of metabolites of the secondary metabolism.

Novel Aspect

Novel method collection for analyzing exudates.

We demonstrate the feasibility of analyzing exudates under field conditions by mass spectrometry and the benefits of Eco Metabolomics study.

References

1. Chaparro J. M., Badri D. V., Bakker M. G., Sugiyama A., Manter D. K., Vivanco J. M., *Journal of Biological Chemistry*, 288, 4502-12 (2013)
2. Lamont B. B., *Plant and Soil*, 2003. 284, 1-19 (2003)
3. Bais, H. P., Park S. W., Weir T. L., Callaway R. M., Vivanco J. M., *Trends in Plant Science*, 9, 26-32 (2004)
4. Bais H. P., Weir T. L., Perry L. G., Gilroy S., Vivanco J. M., *Annual Reviews*, 57, 233-266 (2006)
5. Strehmel N., Bottcher C., Schmidt S., Scheel D., *Phytochemistry*, 30,35-46 (2014)
6. van Dam N. M., Bouwmeester H. J., *Trends in Plant Science*, 21, 256-65 (2016)
7. Weir T. L., Park S. W., Vivanco J. M., *Current Opinions in Plant Biology*, 7, 472-479 (2004)

TP-202 / A NOVEL LABELED METABOLOMICS WORKFLOW APPLYING ISOTOPE RATIO OUTLIER ANALYSIS (IROA) AND SWATH® ACQUISITION FOR UNAMBIGUOUS COMPOUND IDENTIFICATION

Papan Cyrus⁽¹⁾ - *Chris Beecher*⁽²⁾ - *Felice A. de Jong*⁽²⁾ - *Baljit K. Ubhi*⁽³⁾ - *Jörg Dojahn*⁽⁴⁾

SCIEX Germany GmbH, Sciex, Darmstadt (1) - IROA Technologies, University Florida, Bolton MA (2) - SCIEX USA, Sciex, Redwood City (3) - SCIEX Germany (4)

Introduction

Metabolomics focuses on the chemical processes central to cellular metabolism. Mass spectrometry and specifically data dependent workflows tend to be the choice for the measurement of these metabolites. Data independent techniques such as SWATH® Acquisition are different in that they allow for unbiased data collection and MSMS of every single mass precursor can be collected allowing for information rich datasets. However, unambiguous metabolite identification can be increasingly challenging due to the lack of databases, chemical noise and isobaric compounds. The SWATH analysis of the Isotope Ratio Outlier

Analysis (IROA) labeled Internal Standard (IS) provides the first mechanism for simultaneous and unambiguous compound identification and quantitation for unbiased metabolomics analysis.

Methods

A biochemically complex Internal Standard (IS) which contains 100's of biochemicals, each with an IROA isotopic pattern, was added to clinical samples to accurately identify and quantitate complex mixtures without the need for baseline separation, and overcome variances introduced sample-to-sample or by ion suppression. SWATH fragmentation of the IROA peaks completely differentiates fragments, and artifacts. The identification of all IROA compounds and their fragments by ultra-high-resolution mass measurement make it possible to determine the empirical formula for all fragments. Data were collected a TripleTOF® 6600 System in SWATH® Acquisition using a variable window strategy which defined windows of varying mass ranges to be applied in areas of the chromatogram where there are many ions co-eluting.

Preliminary data

Applying a variable window SWATH Acquisition strategy to an IROA Internal Standard (IS) spiked sample made it possible to unambiguously identify and accurately quantify hundreds of biochemicals in a single unbiased metabolomics analysis using a TripleTOF® high resolution mass spectrometer. The IS contains 500+ well characterized compounds, which migrated in an HPLC separation with their natural abundance isotopomers, to provide for both identification and standard quantities for accurate measurement even in a non-baseline separated, "unbiased" metabolomics separation. Using traditional DIA, all compounds with the same retention time are fragmented without selection. SWATH DIA subjects all ions within a discrete m/z window to fragmentation allowing specific precursor ions to be selected, making it easier to analyze fragmentation spectra. However, a corresponding spectral library of metabolites is required for accurate identification. Here, we present SWATH-IROA DIA whereby uniquely-labeled IROA metabolites were captured within discrete SWATH windows, and subjected to fragmentation. IROA fragments and adducts were shown to have the identical labeling patterns of their precursor ions, with defined formulae. All artefactual (non-IROA) peaks from the SWATH window were eliminated and data was quantitated based on MSMS peaks. The combination of IROA and SWATH allow a path in which a basic metabolomic-style system may be used for the accurate clinical quantitation of several hundred compounds in a single sample without the need for a base-line separation. Specific software was developed to automatically find, quantitate and identify all natural-abundance peaks that corresponded to their known IROA isotopomers. The identification of compounds of unknown identity is simplified because all fragments are identified by their complete formula making the mode of fragmentation of the parent compound clear. Fragments are identified as they share the daughter fragments of their parent compounds.

TP-203 / UNTARGETED METABOLOMICS IN PLANT BIOCHEMISTRY

Stefanie Döll - Hendrik Treutler - Steffen Neumann - Dierk Scheel

Leibniz Institute of Plant Biochemistry, Stress and Developmental Biology, Halle

Keywords: Untargeted metabolomics, LC-MS, plant science, substance annotation, functional genomics

Introduction

Search for differentially accumulated metabolites in plant science.

A frequent task not only in plant sciences to find biomarkers for certain conditions ("treatment vs. control"), find chemical differences and similarities between cultivars, subspecies or species and to aid functional genomics by determining the effects of single mutations on the metabolome. This is done by untargeted metabolomics, a technique that works without prior knowledge to the metabolites' identity.

Methods

Experimental setup for UPLC-HR-Q-ToF MS.

For this task, we established a work flow starting from an appropriate experimental setup including numerous quality checks for liquid chromatography coupled to High Resolution Q-ToF MS and MS/MS. This is followed by data analysis with XCMS [1] and MetaboScape (Bruker) including appropriate statistical methods. A combination of tools is applied for substance annotation and identification.

Results

A combination of tools is a start to a nearly comprehensive annotation.

Our approach extracts substances of significance with a very high confidence. Limited by the extraction protocol and LC method, the focus is on semipolar compounds. After determining relevant features, the challenge is the identification of these unknown metabolites: In-house libraries are the basis for batch processing of whole feature tables in MetaboScape. This batch search is extended by publicly available databases (Massbank [2], GNPS [3], ReSpect). Further annotation is done using the in-silico fragmentation tool MetFrag[4]. It is clear that most of the estimated 200,000 metabolites in the plant kingdom are not yet included in databases and therefore, cannot be annotated easily. But, in many cases it is already sufficient to know the compound classes of differentially accumulated metabolites. MetFamily [5] groups features with available MS/MS information according to common fragments or neutral losses and allows a survey of regulated compound classes.

Conclusions

A workflow for automated annotation of known and clustering of unknown metabolites.

Advances in processing large metabolomics data sets, the growing wealth of publicly available databases and the development of customized tools for metabolite annotation largely facilitates the batch annotation of previously known and unknown metabolites of biological relevance.

Novel Aspect

A comprehensive approach to annotate metabolites of biological relevance.

References

1. Smith C. A., Want E. J., O'Maille G., Abagyan R. Siuzdak G., Anal. Chem. 78, 779-787(2006)
2. Horai H. et al., J. Mass Spectrom. 45, 703-714 (2010)
3. Wang M. X. et al., NatBiotechnol. 34, 828-837 (2016)
4. Ruttkies C., Schymanski E. L., Wolf S., Hollender J. Neumann S., J. Cheminform. 8, 3 (2016)
5. Treutler H., Tsugawa H., Porzel A., Gorzolka K., Tissier A., Neumann S., Balcke G. U., Anal. Chem. 88, 8082-8090 (2016)

TP-204 / AN INVESTIGATION OF BIOMARKERS IN THE SERUM OF DEVELOPMENTAL DISORDERS USING AN INTELLECTUAL DISABILITY MODEL MOUSE, WHICH SHOW IMPAIRED SOCIAL BEHAVIOURS

Keiji Gamoh⁽¹⁾ - Mai Nishimura⁽¹⁾ - Takumi Akashi⁽¹⁾ - Yoshinori Nishiwaki⁽¹⁾ - Taiki Kashio⁽²⁾ - Shinichi Mitsui⁽²⁾

Chemistry Department, Faculty of Education, Kochi University, 780-8520, Kochi, Japan (1) - Department of Rehabilitation Sciences, Gunma University Graduate School of Health Sciences, 371-8514, Maebashi, Japan (2)

Keywords: Developmental Disorder Model, Motopsin Knockout (KO) Mice, Biomarkers, Serum, High Resolution LC/MS

Introduction

Our aim is identifying biomarkers of developmental disorders by a liquid chromatography/mass spectrometry (LC/MS). In the present study, serum samples of motopsin KO mice were used for the

exhaustive analysis of biological metabolites-amino acids. We demonstrated the exhaustive analysis of mouse serum using a reversed-phase separation mode and an electrospray ionization method of the LC/MS and compared metabolites between motopsin KO and wild-type (WT) mice.

Methods

Acetonitrile and 1% acetic acid, as well as ^{15}N -glutamic acid as an internal standard, were added to serum samples of motopsin KO and WT mice, and the solution was centrifuged. After the solid-phase extraction of the supernatant using a ODS cartridge, the eluent was injected to the Intra Amino Acid column (Imtakt). The high resolution LC/MS system is composed by Accela (LC) and Orbitrap Exactive (Thermo Fisher Scientific).

Results

Various kinds of amino acids, neurotransmitters, and their metabolites were simultaneously separated and detected by the LC/MS method using standard reference materials. To analyze the effects of genotype and sex, we performed two-way analysis of variance to compare serum concentrations of those substances between KO (male, $n=5$; female, $n=3$) and WT mice (male, $n=4$; female, $n=5$). Serum concentrations of Glu, Gly, and Asn were significantly higher in KO mice than WT mice ($p=0.016$, $p=0.028$, $p=0.010$, respectively), whereas that of DOPAC was lower in KO mice than WT mice ($p=0.003$). Interestingly, post hoc analyses revealed that such difference was significant in males, but not females. Asp, Ser, dopamine, and DOPAC showed higher concentration in males than females ($p=0.006$, $p=0.044$, $p=0.000$, $p=0.000$, respectively). Glu and Gly are neurotransmitters and Asn affects the metabolism control of the cell function of the nervous system and the brain tissue. DOPAC is a metabolite of the neurotransmitter dopamine and is degraded to form HVA. Interestingly, the current result is consistent with our previous research demonstrating significantly higher concentration of Glu in saliva of autism children compared to typically developing children. According to the above results, possible biomarkers of developmental disorders are suggested. It will be necessary to increase the number of subjects for identifying biomarkers with statistical reliability.

This work was supported in part by JSPS KAKENHI Grant No. 17K18650.

Novel Aspect

We demonstrated a possible strategy to find biomarkers of a developmental disorder based on a comparative LC/MS analysis of serum samples between motopsin KO and WT mice.

References

1. S. Mitsui et al., Eur. J. Neurosci. 30, 2368-2378 (2009).
2. K. Suzuki, Shoni no Seishin to Shinkei (in Japanese) 51, 133-142 (2011).
3. C. J. McDougle, C. A. Erickson, K. A. Stigler, J. Clin. Psychiatry 66, 9-18 (2005).
4. S. H. Fatemi, A. R. Halt, J. M. Stary, Biol. Psychiatry 52, 805-810 (2002).

TP-205 / ELECTROCHEMICAL OXIDATION OF PRIMARY BILE ACIDS: SIMULATION OF ENDOGENOUS METABOLISM

Laura Navarro Suarez - Lea Brückner - Sascha Rohn

Hamburg School of Food Science, University of Hamburg, Hamburg

Keywords: Electrochemistry, Oxidation, Bile Acids, EC-MS

Introduction

Electrochemical oxidation is a rapid, relatively inexpensive approach to simulate natural metabolic redox events and potential subsequent reactions in vitro [1,2]. This makes it possible to study biologically relevant

oxidative and reductive processes in an organism and thus, to replace animal studies within the meaning of the directive 2010/63/EU [3]. In this case, coupling of an electrochemical and a mass spectrometric system should simulate phase I and phase II metabolism of bile acids.

Methods

Electrochemical oxidation was performed using a ROXY™ EC system (equipped with a μ -Prep cell consisting of a BDD working electrode, a titanium counter electrode and a Pd/H₂ reference electrode). A solution containing 100 μ M CA/CDCA was injected into the EC-system. First scan mode was applied (0 V to 3 V in 0.3 V steps). Afterwards a constant potential was used. Detection of oxidation products was performed with an ESI-MS ion trap mass analyzer in positive ion mode.

Results

First, the applicability of an electrochemical oxidation for transforming bile acids was evaluated. Simulation of endogenous metabolism of primary bile acids was successful. It was shown that it is possible to oxidize cholic acid and chenodeoxycholic acid under various conditions. Best parameters were acidic pH values and a constant potential (DC-mode). It was also shown that switching to pulse-mode (as mentioned in literature [4]) did not increase the yield. Moreover, it was possible to generate mg-amounts of the oxidized products with a SynthesisCell™ and a general comparison to mass spectra of commercially available bile acids standards.

Conclusions

Electrochemical oxidation using a ROXY™ EC system coupled to a mass spectrometric system is suitable to simulate phase I and phase II metabolism of primary bile acids. It is possible to investigate biosynthetic pathways without using animals or animal tissues.

Novel Aspect

Applicability of an electrochemical oxidation of cholic acid and chenodeoxycholic acid using a ROXY™ EC system was successfully demonstrated.

References

1. Oberacher H., Pitterl F., Chervet J. P., LC GC Europe 28, 138-150 (2015).
2. Lohmann W., Baumann A., Karst U., LC GC Europe 23, 8-16 (2010).
3. Faber H., Vogel M., Karst U., Anal. Chimica Acta 834, 9-21 (2014).
4. Nouri-Nigjeh E., Permentier H. P., Bischoff R., Bruins, A. P., Anal. Chem. 83, 5519-5525 (2011).

TP-206 / LOW INVASIVE METABOLOMIC PHENOTYPING: A POTENTIAL ALTERNATIVE FOR GENOTYPING OF GLIOMAS

Paulina Zofia Goryńska⁽¹⁾ - *Krzysztof Goryński*⁽¹⁾ - *Kamila Chmara*⁽¹⁾ - *Jacek Furtak*⁽²⁾ - *Dariusz Paczkowski*⁽²⁾ - *Marek Harat*⁽²⁾ - *Barbara Bojko*⁽¹⁾

Faculty of Pharmacy Collegium Medicum, Department of Pharmacodynamics and Molecular Pharmacology, Nicolaus Copernicus University in Toruń, Bydgoszcz (1) - 10th Military Research Hospital, Department of Neurosurgery, Bydgoszcz (2)

Keywords: Glioma, SPME, diagnostic method, mass spectrometry

Introduction

Diagnostic of human brain tumors has developed greatly in the recent few years. Proper classification of tumor type plays significant role in prediction of the response to a given treatment, but histology often fails in this regards and genotyping is used instead. However, several mutations can affect therapy efficiency

and identification of molecular biomarkers could be useful complement to current testing.

Methods

In the current study, described approach was applied to human brain tumors. To provide minimal tissue damage ca. 0.2 mm diameters probe had been used to extract metabolites from brain tumor sample. For chromatographic analysis two types of stationary phase had been used: pentafluorophenyl (PFP) and underivatized silica (HILIC). The samples were analyzed on Q-Exactive Focus mass spectrometer and Compound Discoverer 2.1 were used for data processing and statistical analysis.

Results

Metabolomics studies with SPME-LC-HRMS showed that there are several compounds which discriminate tumors of different origin as well as gliomas of different malignancy. Identification of these features revealed involvement of creatine, hypotaurine, cystathionine, etc. in cancerogenesis of gliomas, which was already reported in the literature based on cell line study. It was investigated, that higher concentration of other from those above mentioned naturally produced in the body aminoacids and its homologues indicates the presence of high grade of malignancy glioma in examined patient. Moreover in current study it has been shown that there are several endogenous metabolites that differentiate high grade of malignancy glioma in patients with genetic mutations (IDH1, IDH2, metylation, deletion 1p19q) from those without that mutation. Comparison of the data obtained from genotyping and metabolomics profiling (phenotyping) exhibited very good correspondence.

Conclusions

According to performed experiments in the current study, the use of SPME fibers for extraction and separation via liquid chromatography coupled to high resolution mass spectrometer provided an appropriate workflow for comprehensive brain tumor metabolomics analysis. The simplicity of the approach when compared to routinely used diagnostic methods indicates a great potential of the method to serve as an intra-operative tool in the future.

Novel Aspect

The experimental results obtained using low invasive sample preparation method coupled with LC-MS provide new insight into the diagnostic options of human brain tumors.

TP-207 / $^{13}\text{CO}_2$ ISOTOPE LABELLING OF ARABIDOPSIS THALIANA: TRACKING LABELLED CARBON IN METABOLITES FROM LEAVES TO ROOTS TO EXUDATES

Karin Gorzolka⁽¹⁾ - Gerd Balcke⁽²⁾ - Hendrik Treutler⁽¹⁾ - Steffen Neumann⁽¹⁾ - Dierk Scheel⁽¹⁾

Leibniz Institute of Plant Biochemistry, Stress and Developmental Biology, Halle, Saale (1) - Leibniz Institute of Plant Biochemistry, Cell and Metabolic Biology, Halle, Saale (2)

Keywords: $^{13}\text{CO}_2$ isotope labelling, metabolites, Arabidopsis thaliana tissues, GC-MS, LC-MS

Introduction

How fast is carbon transferred from leaves to roots and exuded?

Plants interact with their belowground environment by the secretion of compounds, a process called exudation [1]. Little is known about this exudation process, exudate biosynthesis and upstream metabolic dynamics. To resolve metabolite origins, we used $^{13}\text{CO}_2$ to track labelled carbon in metabolites from fixation by photosynthesis via transport into roots down to exudation with respect to time dependent dynamics covering light and dark phases.

Methods

^{13}C labelling and spatio-temporal metabolite profiling of Arabidopsis thaliana.

Arabidopsis thaliana plants were pulse labelled in a $^{13}\text{CO}_2$ atmosphere. Aerial plant parts, roots and exudates were collected 3, 6, 9, 20 and 23 hours after labelling and analyzed for primary and secondary metabolites using GC-MS and HPLC-MS. Metabolite labelling was quantified by the investigation of shifts in the isotope pattern in the MS. An XCMS based analysis pipeline was developed for the detection and quantification of isotopes in GC-MS and LC-MS data.

Results

Arabidopsis thaliana - Labelling dynamics of metabolites.

In GC-MS, up to 14000 peaks were detected and clustered in isotopes, resulting in more than 1000 isotope patterns for approximately 350 peak groups (metabolites). Fast labelling was observed in mono- and disaccharides, amino acids and organic acids. Tissue distinction and time-scale harvesting demonstrated that labelling often followed a time dependent delay from leaves (early labelling) to roots to exudates (late labelling); however, spatio-temporal labelling dynamics were still metabolite specific. Especially in saccharides, interesting labelling patterns were observed during the dark phase.

In LC-MS, up to 15 000 peaks were clustered in about 500 peak groups that contained more than 1500 isotope patterns. Labelling was observed in many metabolites, *e.g.* glucosinolates and related compounds in leaves, roots and exudates, providing the opportunity to compare labelling between tissues and metabolites. Co-labelling dynamics pointed to common biosynthetic pathways and metabolite classes and facilitated metabolite identification.

Conclusions

^{13}C tracking in metabolites from leaves to exudates: Biosynthetic connections.

This dataset provides insight into the dynamics between ^{13}C fixation in the aboveground plant, transportation into the roots and exudation of metabolites. The comparison and correlation of labelling patterns revealed biosynthetic connections and consequently facilitated metabolite identification (*e.g.* fatty acid amines). Time dependent sampling resolved changes in metabolite labelling after switches from light to dark phase and back to light.

Novel Aspect

The first MS based study of ^{13}C tracking through leaves, roots and exudates in *Arabidopsis*

We applied MS based ^{13}C labelling analysis of *Arabidopsis* plant tissues and exudates using GC-MS and LC-MS for a comprehensive profiling of primary and secondary metabolites.

References

1. Bais H. P., Weir T. L., Perry L. G., Gilroy S., Vivanco J. M., *Annu. Rev. Plant Biol.* 57, 233 – 266 (2006)

TP-208 / METABOLOMIC APPROACH TO INVESTIGATE METABOLITE ALTERATIONS ATTRIBUTED TO A 25-HYDROXYVITAMIN D IN HEALTHY KOREAN ADULTS

Mi-Ri Gwon⁽¹⁾ - Bo Kyung Kim⁽²⁾ - Boram Ohk⁽²⁾ - Seungil Cho⁽³⁾ - Sook Jin Seong⁽⁴⁾ - Young-Ran Yoon⁽⁵⁾

Cell and Matrix Research Institute, Department of Biomedical Science, School of Medicine, Kyungpook National University, DAEGU (1) - Department of Biomedical Science, School of Medicine, Kyungpook National University, DAEGU (2) - Molecular Diagnostics and Imaging Center, Department of Biomedical Science, School of Medicine, Kyungpook National University, DAEGU (3) - Clinical Trial Center, Kyungpook National University Hospital, DAEGU (4) - Department of Biomedical Science, School of Medicine and Clinical Trial Center, Kyungpook National University and Hospital, DAEGU (5)

Keywords: 25-hydroxyvitamin D [25(OH)D], UHPLC-QTOF/MS, untargeted metabolomics analysis

Introduction

Vitamin D deficiency is prevalent in Korean adults, affecting up to 80% of them [1]. Vitamin D deficiency,

which impairs bone mineralization, may induce many diseases including chronic autoimmune, infectious diseases, cancer, cardiovascular disease, and diabetes. This appears to be associated with levels of the main metabolite, 25(OH)D [2]. We examined the association between serum levels of 25(OH)D and metabolite alterations in healthy Korean adults.

Methods

Forty-five healthy Korean adults were divided into two groups based on serum levels of 25(OH)D (23.0 and 10.4 ng/mL on average for high (n=21) and low (n=24) group, respectively). Demographic, anthropometric, and blood profile parameters were measured. Metabolic profiling was performed using UHPLC-QTOF/MS and multivariate statistical analysis. The candidate metabolites were tentatively identified through the search of web databases.

Results

No significant differences were observed in demographic, anthropometric, and blood profile parameters between two groups. For untargeted metabolomic analysis, 47,499 spectral features were detected from UHPLC-QTOF/MS data using XCMS software. Score plots, obtained from principal component analysis and orthogonal projections to latent structures-discriminant analysis (OPLS-DA), showed a tendency that two groups were separated. R² (goodness of fit) and Q² (predictability) were 0.654 and 0.125, respectively. Metabolites which had the variable importance in projection value ≥ 1.0 and p-value ≥ 0.05 were identified using web databases such as Human Metabolome database (HMDB) and ChemSpider. These results suggest that serum levels of 25(OH)D are related to metabolic alterations in healthy Korean adults. And we will further study to identify putative metabolites and assess their association with synthetic mechanism of 25(OH)D.

Conclusions

This research suggests that serum levels of 25(OH)D have influence on metabolic alterations in healthy Korean adults. However, further research with larger sample sizes and other clinical trial designs (e.g., administration of 25(OH)D to osteoporotic or elderly populations) is required before solid conclusion can be reached.

Novel Aspect

First untargeted metabolomic approach to investigate the influence of serum levels of 25(OH)D on metabolic alterations in healthy Korean adults

References

1. Kim BT, Lee SH, Korean J. Fam. Pract. 3, 6-15 (2013)
2. Gozdzik A, Zhu J, Wong BY, Fu L, Cole DE, Parra EJ, J. Steroid Biochem. Mol. Biol. 127, 405-412 (2011)

TP-209 / AN ANALYTICAL FRAMEWORK FOR THE EXPLORATION OF ENDOGENOUS PEPTIDES INTEGRATING PEPTIDE STRUCTURE PREDICTION AND TANDEM MASS SPECTROMETRY

Chihiro Kawano - Eisuke Hayakawa - Hiroshi Watanabe

Okinawa Institute of Science and Technology, Okinawa

Keywords: Peptidomics, Endogenous Peptide, Tandem Mass Spectrometry, Structure Analysis

Introduction

Peptidomics, an analytical workflow for systematic identification of short, native and endogenous peptides, has been widely used to study peptides in a wide variety of biological matrices. The workflow of peptidomics is largely diverted from proteomics, including LC-MS/MS and spectrum matching algorithm.

However, sensitive and reliable identifications for endogenous peptides are still difficult even with recent mass spectrometry instruments.

Methods

To overcome the difficulties in the conventional peptidomics technologies, we have developed an analytical framework dedicated to the exploration of naturally occurring endogenous peptides. We created a analytical workflow combining tandem mass spectrometry and a series of dedicated informatics modules that enables highly sensitive exploratory search for endogenous peptides.

Results

One of the modules creates potential forms of peptides based on user-defined structural properties such as substructures, motifs, potential cleavage sites, post-translational modification and other molecular properties. The predicted forms of peptides can be used as a target sequence dataset of MS/MS peptide search in a top-down fashion. This approach effectively reduces the search space of MS/MS peptide search, therefore increasing the sensitivity of identification. Moreover, in-silico fragmentation spectra of predicted peptides are generated and their spectral similarities are searched against the experimental fragment spectra acquired from sample. This allows “fuzzy” matching to identify peptides with unexpected post-translational modifications or amino acid substitution with high-sensitivity.

Conclusions

Our analytical framework is built in a modular together with other informatics, data analysis tools, and therefore is flexible and applicable to any kind of peptides in a wide variety of matrices. We discuss the effectiveness and flexibility of our framework by demonstrating identifications of novel endogenous peptides.

Novel Aspect

A novel analytical workflow to improve identification of endogenous peptides by combining tandem mass spectrometry and a series of dedicated informatics tools.

For information please contact: eisuke.hayakawa@oist.jp

TP-210 / TOWARDS THE DECIPHERING OF THE MECHANISM OF ACTION OF AN EXPERIMENTAL DRUG AGAINST CHAGAS DISEASE, USING MULTI-OMICS PLATFORM

Kristin Hennig⁽¹⁾ - *Josephine Abi-Ghanem*⁽¹⁾ - *Andrei Bunescu*⁽¹⁾ - *Xavier Meniche*⁽¹⁾ - *Emeline Biliaut*⁽¹⁾ - *Frédéric Bequet*⁽¹⁾ - *Michael Lewis*⁽²⁾ - *John Kelly*⁽²⁾ - *Stephanie Braillard*⁽³⁾ - *Gilles Courtemanche*⁽¹⁾ - *Eric Chatelain*⁽³⁾

BIOASTER, 40 avenue Tony Garnier, 69007 Lyon, France (1) - Department of Pathogen Molecular Biology, London School of Hygiene and Tropical Medicine, Keppel Street, London WC1E 7HT, UK (2) - Drugs for Neglected Diseases Initiative (DNDi), 15 Chemin Louis-Dunant, 1202 Geneva, Switzerland (3)

Keywords: metabolomics, proteomics, lipidomics, NMR, mass spectrometry

Trypanosoma cruzi (T. cruzi), a kinetoplastid protozoan parasite, is the etiologic agent of Chagas disease. The disease affects millions of people and represents a major burden for the socio-economic development of affected countries (mainly Latin America). The mobility of patients infected with Chagas disease makes it a global public health burden. So far, no vaccine exists and current drugs suffer from significant drawbacks, which limit their use in disease endemic areas. In the foreseeable future, the very few drugs currently in clinical development do not offer new alternatives to patients. Thus, there is a need to identify and develop novel drug candidates with new mechanisms of action to overcome the current situation and eventually feed the Chagas clinical trials pipeline.

The objectives of our study were to identify metabolites and proteins markers of the efficacy of an experimental drug, S205, and to determine the related metabolic pathways to formulate hypotheses on its unknown mode of action. To this aim, we focused on uninfected and *T. cruzi* infected rat myoblast cells, in presence or not of the drug targeting intracellular amastigotes form of the parasite (intracellular phase in the life-cycle of trypanosomes). The drug effect on cells was analyzed after 24h of exposition. We analyzed polar extracellular metabolites by 600 MHz NMR and intracellular lipids and proteins were determined by LC-HRMS using a Q-Exactive (Thermo) and a Q-TOF (Bruker Maxis HD), respectively.

Through our profiling approach, we identified potential markers for infection and drug mode of action. Indeed, the GAPDH (glyceraldehyde-3-phosphate dehydrogenase) protein of *T. cruzi* was identified as a potential target of the drug S205 by proteomics and metabolomics has shown that the drug affects the energy metabolism and hence, through the impact on Acyl-CoA synthesis, the formation of lipids. Thus, S205 seems to act through a specific target and both technologies applied (proteomics and metabolomics) complement each other in the description of the mechanism.

We have shown the relevance of proteomics and metabolomics approaches to contribute to the elucidation of drug mechanism of action as a complement to classical methods. The methods are applicable even to complex systems such as mammal cells infected with parasites.

TP-211 / METABOLOMIC APPROACH FOR THE ESTIMATION OF AGE OF BLOODSTAIN

Hee-Gyoo Kang - Yoo-Jin Lee - Ae Eun Seok - You-Rim Lee - Arum Park - Sora Mun - Hyo-Jin Kim

Department of Senior Healthcare, BK21 Plus Program, Graduate School, Eulji University, Daejeon

Keywords: Bloodstain, Elapsed time, Forensic, HPLC Q-TOF, Metabolite

Introduction

Bloodstain is created by the blood left at the crime scene, and it has a lot of information including genetic information. [1]It is very meaningful because it can identify the time of the attack at the incident. Although various methods have been studied to estimate the age of bloodstain, it has not been applied to forensic field. [2-6] In this study, we identify the metabolites that have changed as time passed to discover markers that we can determine the age of bloodstain.

Methods

The bloodstains were prepared using filter paper over a 21-day period. Metabolites were extracted from bloodstain and analyzed using high-performance liquid chromatography-mass spectrometry (HPLC Q-TOF). We performed multivariate analysis, Venn diagrams and partial least squares discriminant analysis (PLS-DA).

Results

We chose 62 candidate molecular features and analyzed partial least squares discriminant analysis (PLS-DA) that the group have been classified with an accuracy of 75.0%, and R² and Q² were 0.7513 and 0.6998, each.

Ten metabolites were identified matching of METLIN database. Ten of Two metabolites, vitamin D3/cholecalciferol and 1-palmitoyl-2-glutaryl phosphatidylcholine (PGPC), are up-regulated as time passed, whereas ergothioneine is decreased.

Conclusions

We identified several candidate metabolites that are up- or down-regulated over time so that we are able to measure the age of the bloodstain.

Novel Aspect

In this study, we aimed to discover metabolites markers so that we can determine the elapsed time of

bloodstain.

References

1. Butler J. M., Forensic DNA typing: biology, technology, and genetics of STR markers; Academic Press, (2005)
2. Anderson S., Howard B., Hobbs G. R., Bishop C. P., A method for determining the age of a bloodstain, 148, 37-45 (2005)
3. Strasser S., Zink A., Kada G., Hinterdorfer P., Peschel O., Heckl W. M., Nerlich A. G., Thalhammer S., Age determination of blood spots in forensic medicine by force spectroscopy, 170, 8-14 (2007)
4. Bremmer R. H., Nadort A., van Leeuwen T. G., van Gemert M. J., Aalders M. C., Age estimation of blood stains by hemoglobin derivative determination using reflectance spectroscopy, 206, 166-171 (2011)
5. Inoue H., Takabe F., Iwasa M., Maeno Y., Identification of fetal hemoglobin and simultaneous estimation of bloodstain age by high-performance liquid chromatography, 104, 127-131 (1991)
6. Li B., Beveridge P., O'Hare W. T., Islam M., The estimation of the age of a blood stain using reflectance spectroscopy with a microspectrophotometer, spectral pre-processing and linear discriminant analysis, 212, 198-204 (2011)

TP-212 / METABOLOMIC ANALYSIS OF HUMAN PLASMA TO CHARACTERIZE METABOLITE CHANGES ASSOCIATED WITH EXPOSURE TO PERSISTENT ORGANIC POLLUTANTS

Bo Kyung Kim⁽¹⁾ - ***Mi-Ri Gwon***⁽²⁾ - ***Min-Hee Kwon***⁽¹⁾ - ***Boram Ohk***⁽¹⁾ - ***Seungil Cho***⁽³⁾ - ***Sook Jin Seong***⁽⁴⁾ - ***Duk-Hee Lee***⁽⁵⁾ - ***Young-Ran Yoon***⁽⁶⁾

Department of Biomedical Science, School of Medicine, Kyungpook National University, Daegu (1) - Cell and Matrix Research Institute, Kyungpook National University, Daegu (2) - Molecular Diagnostics & Imaging Center, Department of Biomedical Science, School of Medicine, Kyungpook National University, Daegu (3) - Clinical Trial Center, Kyungpook National University Hospital, Daegu (4) - Department of Preventive Medicine, School of Medicine, Kyungpook National University, Daegu (5) - Department of Biomedical Science, School of Medicine, Clinical Trial Center, Kyungpook National University, Kyungpook National University Hospital, Daegu (6)

Keywords: persistent organic pollutants, organochlorine pesticides, polychlorinated biphenyls, metabolic syndrome, UHPLC-QTOF/MS, untargeted metabolomics

Introduction

Humans are widely exposed to persistent organic pollutants (POPs) such as organochlorine pesticides (OCPs) and polychlorinated biphenyls (PCBs). Because of their high lipophilicity, POPs are known to bioaccumulate mainly in adipose tissue and to be resistant to biodegradation [1]. Low-dose exposure to OCPs and PCBs is related to metabolic syndrome [2]. Here, we investigated alterations of endogenous metabolites in human plasma caused by OCPs and PCBs exposure.

Methods

A total of 203 subjects (age≥60) participated in this research and serum POPs levels were measured. All subjects were categorized into quartile based on serum concentration levels of OCPs or PCBs. Plasma samples were analyzed using ultra-high performance liquid chromatography (UHPLC)-quadrupole time-of-flight (QTOF)/mass spectrometry (MS), then multivariate statistical analysis was performed. The metabolite features of variable importance in projection value (VIP)≥ 1.0 were compared with web databases.

Results

A total of 37 POPs, 4 OCPs and 6 PCBs were selected and considered in this study. Through the untargeted metabolic profiling, the 32,301 ion peaks were detected. Score plots obtained from principal component

analysis (PCA) and partial least squares-discriminant analysis (PLS-DA) accounted for the relationship among variables, especially Q1 and Q4. The PLS-DA results showed a tendency to divide that Q1 group is separated from Q4 groups by OCPs concentration, not by PCB's. The 4,491 ion peaks were acquired (VIP \geq 1.0), and considered to contribute separation of two groups. Thirty three compounds were obtained which had a possibility to consider the accurate mass fragmentation pattern. These results suggest that serum OCPs levels can cause metabolic alterations in human. And we will further study to identify putative metabolites that contribute to metabolic changes caused by OCPs exposure.

Conclusions

Through this study, we will show UHPLC-QTOF/MS-based metabolomics research can be a useful tool to evaluate the effects of POPs exposure by characterizing the metabolic alterations in general population by serum OCPs levels and figuring out the relation of POPs exposure and metabolic syndrome.

Novel Aspect

The large scale untargeted metabolomics approach with human plasma to characterize metabolic alterations associated POPs exposure and metabolic disease.

References

1. El-Shahawi MS, Hamza A, Bashammakh AS, Al-Saggaf WT, Taranla, 80, 11587-1597 (2010).
2. Lee DH, Steffes MW, Sjödin A, Jones RS, Needham LL, Jacobs DR Jr., Environ Health Perspect, 118, 1235-1242 (2010).

TP-213 / A RAPID MICROBORE METABOLIC PROFILING (RAMMP) ANALYTICAL PLATFORM FOR DISCOVERY METABOLOMICS AND LIPIDOMICS

Adam King⁽¹⁾ - **Lauren Mullin**⁽¹⁾ - **Ian Wilson**⁽²⁾ - **Paul Rainville**⁽¹⁾ - **Robert Plumb**⁽¹⁾ - **Lee Gethings**⁽¹⁾ - **David Heywood**⁽³⁾ - **Garth Maker**⁽⁴⁾ - **Robert Trengove**⁽⁵⁾

Waters Corporation, Scientific Operations, Wilmslow (1) - Imperial College, Computational and systems medicine, London (2) - Waters Corporation, Omics business development, Wilmslow (3) - Murdoch University, Medical and Molecular sciences, Perth (4) - Murdoch University, Separations science and Metabolomics laboratory, Perth (5)

Keywords: Metabolomics, Lipidomics, Rapid Phenotyping, UPLC-QToF.

Introduction

Metabolic phenotyping has been employed in large epidemiological studies in the effort to discovery more about human health and disease [1 & 2]. UPLC and hybrid mass spectrometers have been essential tools in analyzing the matrices chosen for these studies [3]. However, long acquisition times per sample (15 – 30 mins) have meant large batches can take days to analyse [2]. Here we have developed a suite of rapid profiling methods reducing acquisition time to ~3 mins per sample.

Methods

Rodent urine and human plasma extracts (1:4 v/v IPA) were analysed using Waters Acquity I-class UPLC system coupled to a Synapt G2-Si QToF capable of ion mobility separation. Urinary metabolites were chromatographically resolved through HILIC separation with lipids by reversed-phase, both utilizing 5.0 cm X 1 mm i.d. columns. Mobile phases consisted of water, acetonitrile for HILIC and the addition of IPA for lipid, both with ammonium formate.

Results

Standard LC-MS methods for metabolite and lipid profiling using 2.1 mm i.d. columns were geometrically

scaled down, reducing column i.d., cutting analysis times, injection volumes and mobile phase flow rate ultimately increasing the linear velocity. Each method scaling resulted in a reduced overall acquisition time to approximately 1/3 of the original method. Compound standard mixtures were used to confirm that chromatographic performance and analyte retention was not compromised with the method scaling. The relative retention time of 3 marker compounds between the standard HILIC separation and the rapid method were comparable, demonstrating a preservation of the HILIC mechanism. When analyzing biological matrices, the potential peak capacity of the assay was 28, however, enabling the ion mobility function, this increased to 51 with the overall number of compound ions detected doubled. The inclusion of ion mobility separation and resulting CCS values enabled an additional confirmatory database search for compound identification.

Conclusions

A set of rapid LC-MS methods have been successfully employed in profiling polar and non-polar molecules in rodent urine and human plasma from breast cancer patients. A reduction in analytical run time has enabled large batch data acquisition to be reduced by 60%, eliminating the need for multiple sample batches and the risk of inconsistent results. The inclusion of ion mobility has provided an additional dimension to the data and improved compound identification.

Novel Aspect

Rapid 3 minute long discovery phenotyping assays in conjunction with high resolution accurate mass and ion mobility spectroscopy.

References

- 1 Holmes E., Wilson I D., Nicholson J K., Cell 134, 714-717 (2008).
2. Lewis L R., Pearce J T M., Spagou K., Green M., Dona A C., Yuen A H Y., David M., Berry D J., Chappell K., Horneffer-van der Sluis V., Shaw R., Lovestone S., Elliott P., Shockcor J., Lindon J C., Cloarec O., Takats Z., Holmes E., Nicholson J K., Anal. Chem. 88, 9004-9013 (2016).
3. Rainville P D., Wilson I D., Nicholson J K., Isaac G., Mullin L., Langridge J I., Plumb R S., Anal. Chim. Acta 982, 1-8 (2017).

TP-214 / RAPID AND SIMULTANEOUS DETERMINATION OF 36 AMINO ACIDS IN NON-INVASIVE AND INVASIVE BIOLOGICAL SAMPLES FROM OBESITY PATIENTS BY LIQUID CHROMATOGRAPHY-TANDEM MASS SPECTROMETRY AFTER DERIVATIZATION

Lucyna Konieczna⁽¹⁾ - ***Anna Krawczyńska***⁽¹⁾ - ***Janusz Siebert***⁽²⁾ - ***Piotr Gutknecht***⁽²⁾ - ***Maria Skrzypkowska***⁽³⁾ - ***Łukasz Kaska***⁽⁴⁾ - ***Justyna Bigda***⁽⁴⁾ - ***Maria Proczko-Stepaniak***⁽⁴⁾ - ***Tomasz Bączek***⁽¹⁾

Department of Pharmaceutical Chemistry, Medical University of Gdańsk, Poland (1) - Department of Family Medicine, Medical University of Gdańsk, Poland (2) - Department of Immunology, Medical University of Gdańsk, Poland (3) - Department of General, Endocrine and Transplant Surgery, Medical University of Gdańsk, Poland (4)

Keywords: exhaled breath condensate (EBC) samples, amino acids, derivatization step, LC-MS/MS, bariatric treatment

Introduction

The analysis of free amino acids (AAs) in biological samples is important for diagnosing the health of individuals, because their concentrations are known to vary with various diseases [1-4]. In this study, reliable and high-throughput analytical method for AAs analysis in different biological matrices was proposed and applied to assess their concentration in noninvasively collected breath condensate samples obtained from patients before bariatric treatment.

Methods

A robust and sensitive method for identification and quantification of free amino acids in exhaled breath condensate, plasma, saliva and urine samples, based on liquid chromatography-electrospray ionization-tandem mass spectrometry (LC-ESI-MS/MS) has been developed. Due to the presence of highly polar groups (-NH₂, -OH, -COOH) and their non-volatile nature, amino acids need to be derivatized with propyl chloroformate before LC-MS/MS analysis.

Results

To evaluate the efficiency of the present method, the method was applied to human exhaled breath condensate, urine, plasma and saliva samples collected from obesity patients qualified for surgical bariatric treatment. Adequate preoperative preparation and postoperative care of patients requires profiling of amino acids to control of diet and search potential biomarker of obesity-related diseases. In our developed method, the separation of 36 of physiological amino acids was achieved in 18 min. It was observed that arginine, glutamine, serine, proline, hydroxyproline, glycine, threonine, alanine, proline, methionine, aspartic acid, histidine, lysine, norvaline, glutamic acid, tryptophan, leucine, phenylalanine, isoleucine, asparagine, cystine, cysteine and tyrosine levels are higher in plasma and saliva with reference to breath condensate samples. Only valine concentration are comparable in both saliva and EBC samples.

Conclusions

Propyl chloroformate derivatization followed by LC/MS/MS analysis provides an effective platform for quantitative analysis of amino acids in complex biological samples such as plasma, urine, saliva and exhaled breath condensate samples. The proposed method is rapid, sensitive and can easily be performed in aqueous media. Further it does not require any pre-purification, cleanup and/or lyophilization steps before analysis. The method developed has wide applications for the routine analysis of amino acids in non-invasive breath samples, without the need for blood collection.

Novel Aspect

A sensitive and direct analytical method has been proposed for the first time for simultaneous determination of 36 amino acids in non-invasive EBC samples.

References

1. Chace D. H., Hillman S. L., Lorey F., Cunningham G. C., Clin. Chem. 44, 2405-2412 (1998).
2. Wang T. J., Larson M. G., Vasan R. S. et al., Nat. Med. 17, 448-454 (2011).
3. Yosida H., Kondo K., Yamamoto H. et al., J. Chromatogr. B 998-999, 88-95 (2015).
4. Konieczna L., Pyszka M., Okońska M., Niedźwiecki M., Bączek T., J. Chromatogr. A1542, 72-81 (2018).

TP-215 / IDENTIFICATION OF INTERNAL STANDARD METABOLITES FOR STANDARDIZATION BLOODSTAINS

Jiyeong Lee⁽¹⁾ - **Hyo-Jin Kim**⁽²⁾ - **You-Rim Lee**⁽²⁾ - **Ae Eun Seok**⁽²⁾ - **Arum Park**⁽²⁾ - **Sora Mun**⁽²⁾ - **Yoo-Jin Lee**⁽²⁾ - **Hee-Gyoo Kang**⁽¹⁾

Department of Biomedical Laboratory Science, College of Health Sciences, Eulji University, Seongnam-si, Gyeonggi-do (1) - Laboratory of Signal Transduction and Disease Biomarker Discovery, Department of Senior Healthcare, BK21 Plus Program, Graduate School, Eulji University, Daejeon (2)

Keywords: Bloodstain, Forensic science, Internal standard, Metabolite

Introduction

Bloodstains obtained from crime site provide important information for scene reconstruction. Various researches have been conducted to obtain information of the crime by analyzing bloodstain components. However, the components vary with bloodstain volume and other unpredictable factors. In this study, we

identified internal standard metabolites that having stable volume even long after its deposition for standardize bloodstains.

Methods

Blood was collected from six subjects. Blood stains were kept until following time points: day 0, day 7, day 14, day 21 and day 28. The metabolite profiling was conducted by using a liquid chromatography-tandem mass spectrometry (LC-MS/MS). Coefficient of variation(CV) of molecular features was calculated for each criterion: time point, subject, all data. Fold change and statistical analysis were performed to select suitable candidates for internal standard.

Results

Five molecular features with average of $CVs \leq 5\%$ were selected as candidates for internal standard. Partial least squares discriminant analysis (PLS-DA) and principal component analysis(PCA) shows well grouping of time condition. This results indicated the candidates unrelated to elapsed time. Fold change value of those candidates' abundances was calculated to each time points versus day 0. The stigmaterol had the most stable pattern and the L-methionine remained stable until day 14 and after day 21. Therefore, two metabolites of candidates were suitable for internal standard.

Conclusions

In this study, we selected the candidates of internal standard to obtain the information of initial volume of bloodstain for more accurate analysis. Further studies should be conducted to validate the candidates using more samples and conditions associated with the real field situation. This study provides the basic data to obtain absolute quantitative information and proposed the possibility of bloodstain standardization method using the internal standard.

Novel Aspect

The internal standard metabolites will contribute to establish standardization of bloodstains and help to perform accurate bloodstain analysis in the field.

References

1. Bevel, Tom, and Ross M. Gardner. Bloodstain pattern analysis with an introduction to crime scene reconstruction. CRC press, 2008.
2. Hara, M., et al. "Effects of storage method on DNA degradation in old bloodstain samples. " Forensic Science International: Genetics Supplement Series 5 (2015), e39-e41.
3. Huang Y, Yan J, Hou J, Fu X, Li L, Hou Y (2015) Developing a DNA methylation assay for human age prediction in blood and bloodstain. Forensic Science International: Genetics 17, 129-136
4. Virkler K, Lednev IK (2009) Analysis of body fluids for forensic purposes: from laboratory testing to non-destructive rapid confirmatory identification at a crime scene. Forensic Science International 188 (1-3), 1-17
5. Edelman G, van Leeuwen TG, Aalders MC (2012) Hyperspectral imaging for the age estimation of blood stains at the crime scene. Forensic science international 223 (1-3), 72-77
6. Hanson EK, Ballantyne J (2010) A blue spectral shift of the hemoglobin soret band correlates with the age (time since deposition) of dried bloodstains. PLoS One 5 (9), e12830
7. Bremmer RH, De Bruin DM, De Joode M, Buma WJ, Van Leeuwen TG, Aalders MC (2011) Biphasic oxidation of oxy-hemoglobin in bloodstains. PloS one 6 (7), e21845
8. Bremmer RH, Nadort A, Van Leeuwen TG, Van Gemert MJ, Aalders MC (2011) Age estimation of blood stains by hemoglobin derivative determination using reflectance spectroscopy. Forensic science international 206 (1-3), 166-171
9. Alshehhi S, McCallum NA, Haddrill PR (2017) Quantification of RNA degradation of blood-specific markers to indicate the age of bloodstains. Forensic Science International: Genetics Supplement Series 6, e453-e455

TP-216 / AN INFUSION “SHOTGUN” APPROACH FOR HIGH-THROUGHPUT UNTARGETED METABOLOMICS

Mariateresa Maldini⁽¹⁾ - *Cyrus Papan*⁽¹⁾ - *Zuzana Demianova*⁽¹⁾ - *Jörg Dojahn*⁽¹⁾ – *Baljit K. Ubhi*⁽²⁾

SCIEX, Darmstadt (1) - SCIEX, SCIEX, Redwood City (2)

Keywords: MSMS all, metabolomics, high resolution, DIA

Introduction

Metabolomics studies the complexity and large variety of various chemical compounds. Mass spectrometry-based untargeted metabolomics requires approaches to collect data in an unbiased fashion. Current LC approaches are biased to the column chemistry and metabolites which can be retained on such chemistries. The main goal of researchers in the untargeted metabolomics field is to analyze a large number of samples and obtain the most information in the shortest time with limited sample preparation. Here we evaluate the direct infusion or “shotgun” approach to such analyses by using a technique, Infusion MS/MSALL Acquisition, for untargeted metabolomics. This data-independent technique allows the MS/MS of all possible candidates, improve identifications significantly, allows to have quantitative data and enables retrospective analysis of the data.

Methods

Data were acquired using a TripleTOF® 6600 System (SCIEX) coupled to a high flow liquid chromatography system (Shimadzu) and employing the Infusion MSMSALL Acquisition to collect one survey scan and then MS/MS of every single precursor within a given mass range. Mass ranges were evaluated from 50 to 900 in varied sizes to optimize the number of MS/MS and time needed to collect the data. Data were collected on a small case study of urine extract and statistical analyses was performed.

Data were processed using MasterView™ Software and the accurate mass metabolomics spectral library for fast metabolite identification and confirmation. Statistical analysis was completed in MarkerView™ Software and any quantitation results were generated by MultiQuant™ Software.

Results

By optimizing for mass range, MSMS ALL method was developed to capture all relevant metabolites per ionization mode (positive and negative mode). Metabolites are extracted from several biological samples and directly infused by using the developed method. The Total Ion chromatogram shows the presence of a large number of separate peaks. The extracted ion chromatograms show that most peaks represent, as expected, different metabolite classes, which are characteristic for each sample. Preliminary, by using the fragment filter tool, metabolites classes occurring in the different samples, were identified on the basis of their diagnostic fragment peaks and/or neutral loss. Statistical analysis is successively performed to differentiate samples and highlight rapidly marker compound in each sample. Using the high quality of MS/MS spectra thanks to SWATH acquisition, identification of marker and occurring metabolites is obtained together with the relative amount.

“Shotgun” approach for Metabolomics results a powerful and rapid tool to obtain huge and complementary data of numerous samples even if metabolites present only at very low levels in biological samples.

Post-acquisition data processing software, allows the user to have rapidly structural information and/or identification and quantification of metabolites occurring in complex samples to derive important biological conclusions.

Conclusions

By using Infusion MSMSALL with SWATH acquisition, a rapid metabolic fingerprint of each biological sample can be captured. The optimized method is able to capture all relevant metabolites per ionization mode (positive and negative mode). Verifying the method through a simple biological extract case study highlighted the metabolites responsible for differentiating the different sample groups.

Novel Aspect

An Infusion “Shotgun” Approach for High Throughput Untargeted Metabolomics

TP-217 / SWATH® ACQUISITION IMPROVES METABOLITE COVERAGE OVER TRADITIONAL DATA DEPENDENT TECHNIQUES FOR UNTARGETED METABOLOMICS

Mariateresa Maldini⁽¹⁾ - *Zuzana Demianova*⁽¹⁾ - *Cyrus Papan*⁽¹⁾ - *Jörg Dojahn*⁽¹⁾ - *Baljit K. Ubhi*⁽²⁾

SCIEX, SCIEX, Darmstadt (1) - SCIEX, SCIEX, Redwood City (2)

Keywords: Untargeted metabolomics, SWATH®, Data Dependent Acquisition

Introduction

Data independent acquisition (DIA) workflows are well adopted in quantitative discovery proteomics, but still not commonly used in discovery metabolomics. Data dependent acquisition (DDA) techniques are heavily employed in the field of metabolomics and workflows on mass spectrometers have been adapted so that as much data as possible can be captured. Researchers were limited by the speed of their QTOF mass spectrometers meaning a multiple injection workflow. Also, the stochastic nature of data dependent workflows often means MSMS of low abundant metabolites are often missed. Here, we describe how DIA enables the identification of a higher number of metabolites for untargeted metabolomics workflows compared to traditional DDA approaches thus enabling a broader profile of the metabolome.

Methods

Urine and plasma were processed according to standard extraction protocols. Urine was diluted with water at a ratio of 1:4 (v/v) and centrifuged for prior analysis, whilst plasma was extracted 1:4 (v/v) with ice-cold methanol allowing for protein precipitation. Separation was performed on reverse phase chromatography. DDA and DIA specific settings were chosen and evaluated. The data were acquired on a QqTOF mass analyzer. For the DDA acquisition, we selected the top 5, 10, 15, 20 and 25 precursor ions for MSMS. For DIA, we applied 15, 20 and 30 mass windows with either fixed window (fw) or variable window (vw) widths. results were evaluated by the highest number of identifications and coverage of metabolites in plasma and urine extracts.

Results

At the DDA level, the data demonstrate a significant improvement of metabolite coverage at the MSMS level when comparing the top5 to the top25 DDA method. We show over 100% increase of metabolite coverage in plasma extracts by increasing the number of selected precursor ions for DDA acquisition from top5 to top25. This result highlights the capability of the QqTOF mass analyzer for fast MSMS acquisition, which allows for the fragmentation of a large number of precursors in a single DDA cycle, leading to a larger number of metabolites identified.

In the second part of this study, we evaluated the DIA strategy with various fixed (fw) and variable window (vw) sizes with similar cycle time in a plasma extract. Increasing the number of fixed windows resulted in ~30% gain in metabolite coverage. Using the variable window method resulted in a ~70% gain in metabolite coverage.

We finally applied these experimental approaches to common matrices used in metabolomics studies, namely urine and extracted plasma. We show that a DIA approach applying 20 variable windows can identify up to 55% more metabolites than a traditional top20 DDA acquisition (in a urine extract). More confident MSMS based identifications lead to higher quantifiable metabolites in a metabolite expression experiment, which at the end allows better understanding of the biology. When comparing the performance in extracted plasma it can be observed that applying a DIA approach with 20 variable windows allows significant gains in metabolite coverage (around 55%) versus the top20 DDA acquisition, similar gains as seen in the urine extract.

Conclusions

Data Independent Acquisition allows to quickly obtain a comprehensive identification and quantification of the metabolites in your sample since is able to acquire full scan (MS1) and MS/MS of every single metabolite in a single injection. Variable Window Calculator can be used to optimize Q1 isolation window pattern to achieve the right balance of metabolite coverage and specificity.

Novel Aspect

Data independent acquisition improves metabolite coverage over traditional data dependent techniques for untargeted metabolomics.

TP-218 / ANALYTICAL ION SPECTROSCOPY - A NEW ROUTE TO METABOLITE IDENTIFICATION

Jonathan Martens⁽¹⁾ - Giel Berden⁽¹⁾ - Rianne van Outersterp⁽¹⁾ - Leo Kluijtmans⁽²⁾ - Ron Wevers⁽²⁾ - Jos Oomens⁽¹⁾

FELIX Laboratory, Radboud University, Nijmegen (1) - Translational Metabolic Laboratory, Radboud University Medical Centre, Nijmegen (2)

Keywords: Infrared ion spectroscopy, metabolite identification

Introduction

Small molecule identification is a key component of many areas of the(bio)analytical sciences, including metabolomics and drug development. Here, we demonstrate how infrared (IR) ion spectroscopy is being used for small molecule identification by metabolomics labs. We present technological developments that enable separation by (ultra) high pressure liquid chromatography followed by structural identification by IR spectroscopy of the mass-selected ions (1).

Methods

The experimental setup is based on a modified commercial ion trap mass spectrometer in combination with the tunable infrared radiation of the FELIX infrared free electron laser. Resonant absorption of multiple infrared photons leads to dissociation. Relating the extent of dissociation as a function of the laser frequency gives the infrared spectrum of any mass-selected ion that can finally be matched to reference (or calculated) IR spectra.

Results

The majority of clinical metabolomics laboratories rely on mass spectrometry to identify the presence of abnormal levels of specific small molecules in common body fluids that correlate with diseases, drug treatments, and environmental factors. In a first example, we focus on the differentiation of *N*-acetylmannosamine, a recently discovered biomarker for NANS-deficiency (a new inborn error of metabolism in sialic acid metabolism, from other possible *N*-acetylhexosamines found at the same *m/z* (2). This differentiation is not possible using standard LC-MS protocols available in most metabolomics laboratories. Here, we demonstrate the use of IR ion spectroscopy to cleanly distinguish these three *N*-acetylhexosamines directly from urine and cerebrospinal fluid (CSF). Finally, we show how IR ion spectroscopy has been used to partially identify a previously detected but yet unidentified biomarker of antiquitin deficiency. This peak "X" is thought to be associated with intellectual disability that persists despite current treatment strategies.

Conclusions

The results presented here show how IR ion spectroscopy compares very favorably against alternative techniques, such as NMR spectroscopy, for small molecule identification (3). It maintains the sensitivity (low nanomolar) and (mass) selectivity inherent to mass spectrometry while providing orthogonal

identification on the basis of retention time, mass-to-charge ratio, MS/MS fragmentation patterns and the IR spectral fingerprint.

Novel Aspect

Metabolite identification using combined liquid chromatography-infrared ion spectroscopy

References

1. Martens, J. ; Koppen, V. ; Berden, G. et al., Combined Liquid Chromatography-Infrared Ion Spectroscopy for Identification of Regioisomeric Drug Metabolites. *Anal. Chem.* 89, 4359-4362 (2018).
2. Martens, J. ; Berden, G. ; van Outersterp, et al., Molecular identification in metabolomics using infrared ion spectroscopy. *Scientific Reports*7, 3363, (2017).
3. Martens, J. ; Berden, G. ; Bentlage, H. et al., Unraveling the unknown areas of the human metabolome: the role of infrared ion spectroscopy. *J. Inherited Metab. Dis.*, article online (2018).

TP-219 / MICROFLOW METABOLOMICS DIFFERENTIATES PRE-CLASSIFIED HEALTHY AND CANCER SAMPLES

Dietrich Merkel⁽¹⁾ - *Cyrus Papan*⁽¹⁾ - *Zuzana Demianova*⁽¹⁾ - *Jörg Dojahn*⁽¹⁾ - *Mariateresa Maldini*⁽²⁾ - *Khatereh Motamedchaboki*⁽³⁾ - *Baljit K. Ubhi*⁽⁴⁾ - *Erika Lin*⁽⁵⁾

SCIEX, Germany, Support, Darmstadt (1) - SCIEX, Italy, Support, Milano (2) - SCIEX, USA, Technical Marketing, Redwood City (3) - SCIEX, USA, business unit, Redwood City (4) - SCIEX, USA, Product Management, Redwood City (5)

Keywords: microflow, HILIC, MRM, metabolomics, cancer

Introduction

Identification of metabolites from biological matrices is necessary for studying potential disease biomarkers and to further understand their biological implications. Liquid chromatography–tandem mass spectrometry (LC-MS/MS) analysis is an essential tool for identification and quantitation of metabolites in complex sample matrices. Here, we present HILIC method using a μ flow LC-MRM to evaluate the assay performance on a colon cancer study.

Methods

Metabolite extracts from urine, plasma and an MDCK cell line samples were analyzed with both analytical and μ flow rate on SCIEX ExionLC™ AD HPLC and M3 MicroLC system. HILIC separation was performed using Luna- NH 2 columns (Phenomenex). MRM analysis with pos/neg polarity switching was completed on a QTRAP® 6500+ system. A total of 363 positive and negative ion mode MRM's for a total 312 unique polar metabolites were monitored.

Results

While μ flow LC in general has become increasingly popular, it has not been readily employed in metabolomics because the typically used sample solvent does not allow for injecting larger volumes of samples without sacrificing HILIC chromatographic resolution. However, by simply reconstituting the sample in an organic solvent (95% acetonitrile, pH 9), we could inject up to 5 μ L sample, while maintaining good peak shape. Using this single injection workflow we identified 263 unique metabolites out of 312 from all three matrices, covering all major metabolic pathways. This μ flow HILIC-MRM method on average provided 10 times more sensitivity compared to similar high flow methods. The μ flow HILIC LC-MRM method improved the S/N ratio by up to 60 times for some polar metabolites such as uridine 5'-diphosphate, hippuric acid and cytidine. Excellent chromatographic separation of improved sensitivity and resulted in higher number of detections with a 50% improvement in MDCK cell line, 35% in urine and 11%

in plasma when compared to high flow methods.

Conclusions

For the colon cancer study the previously classified healthy and diseased (Cancer) samples were clearly differentiated by data acquired using this microflow method and data analysis with the multivariate data analysis techniques in MarkerView™ Software. This method identifies the metabolites responsible for the differentiation of the sample grouping (healthy vs cancer), which is clearly being driven by 14 metabolites profile changing due to these conditions.

Novel Aspect

A single microflow targeted metabolomics method allowing detection of over 300 polar metabolites across multiple biochemical pathways.

TP-220 / HIGH-THROUGHPUT LC-MS/MS-BASED CHIRAL METABOLIC PROFILING FOCUSING ON AMINO ACIDS AND RELATED METABOLITES

Yosuke Nakano - Moyu Taniguchi - Eiichiro Fukusaki

Osaka university, Advance Science and Biotechnology, Osaka

Keywords: d-amino acid, chiral metabolomics, enantioseparation, LC-MS/MS

Introduction

Metabolomics has been an evolving science with a wide range of applications in various fields. However, previous studies have rarely focused on metabolite chirality. We previously developed analytical methods for the simultaneous analysis of eighteen chiral proteinogenic amino acids[1, 2]. In this study, to achieve metabolic profiling of chiral amino acids and related metabolites, we developed a high-throughput method using LC-MS/MS.

Methods

LC-MS/MS analysis was performed using the Nexera HPLC System (Shimadzu Corporation, Kyoto, Japan) connected to an LCMS-8060 (Shimadzu) with dual ion source of electrospray ionization and atmospheric pressure chemical ionization in MRM mode. Chromatographic separation was achieved with CROWNPAK CR-I(+)/CR-I(-) (3.0 mm i.d., 150 mm, 5 µm) and CHIRALPAK ZWIX(-) (4.0 mm i.d., 150 mm, 3 µm) (Daicel CPI, Osaka, Japan).

Results

The combination of two types of chiral columns (with binaphthyl-based crown ether and cinchona alkaloid-derived zwitterionic stationary phases) enabled the analysis of 115 chiral and non-chiral metabolites. By finely optimizing MS/MS parameters, the method allowed the highly sensitive (0.001-50 nmol/mL) and wide dynamic range detection (100-50,000) of target analytes in a standard solution without derivatization. We applied the method to food samples (cheese), and successfully quantified trace levels of metabolites such as d-amino acids in samples. Additionally, we performed principal component analysis (PCA) on the metabolome data and obtained unique profiles that reflected metabolite chirality.

Conclusions

This method highlights the highly sensitive detection of trace amounts of targeted metabolites without the need for derivatization. Moreover, we demonstrated the usefulness of our system by simultaneously analyzing both chiral and non-chiral metabolites in samples. According to the profiling, some d- and l-amino acids showed different behaviors in the PCA plot, which could not be observed by conventional methods that rarely consider separation of enantiomers.

Novel Aspect

We anticipate that this method will enable comprehensive profiling to provide novel insights into the discovery of molecular markers utilized in various fields.

References

1. Konya, Y. ; Bamba, T. ; Fukusaki, E., J. Biosci. Bioeng.,121, 349–353 (2016).
2. Nakano, Y. ; Konya, Y. ; Taniguchi, M. ; Fukusaki, E., J. Biosci. Bioeng.,123, 134–138 (2017).

TP-221 / STUDY OF AUXIN BIOSYNTHESIS USING IN VIVO LABELLING

Ales Pencik⁽¹⁾ - *Petr Vonka*⁽¹⁾ - *Martin Kubeš*⁽²⁾ - *Karin Ljung*⁽³⁾ - *Ondřej Novák*⁽¹⁾

Institute of Experimental Botany AS CR & Palacky University, Laboratory of Growth Regulators, Olomouc (1) - Centre of the Region Haná for Biotechnological and Agricultural Research, Palacky University, Department of Chemical Biology and Genetics, Olomouc (2) - Swedish University of Agricultural Sciences, Department of Forest Genetics and Plant Physiology, Umea (3)

Keywords: tandem mass spectrometry, in vivo labelling, auxin, metabolism, Arabidopsis thaliana

Introduction

The phytohormone auxin (indole-3-acetic acid; IAA) plays a fundamental role in various processes of plant growth and development. Crucial for auxin action is the interplay between its transport, biosynthesis, degradation and conjugation. The use of heavy isotope labelling of intermediates in conjunction with mass spectrometry analysis has proven to be a powerful tool for monitoring fluxes through IAA metabolic pathways.

Methods

7-day-old seedlings of Arabidopsis thaliana were incubated with 10 μM [¹³C]indole. Seedlings were harvested after 1, 2, 4, 6 and 24 hours, then extracted and subsequently purified by solid phase extraction (SPE) as described by Novák et al. [1]. Incorporation of [¹³C]indole and its metabolisation in planta were monitored by LC-MS/MS with multiple reaction monitoring transitions and quantified by isotopic dilution method with appropriate internal standards [2].

Results

We performed a feeding experiment using ¹³C labelled indole in Arabidopsis wild type (Col-0) and several IAA biosynthetic and metabolic mutants to investigate auxin metabolism. Incorporation of the ¹³C from labelled precursor into most abundant IAA biosynthetic intermediates (tryptophan, indole-3-pyruvic acid (IPyA), indole-3-acetamide (IAM) and indole-3-acetonitrile (IAN)) as well as metabolites (IAA-aspartate (IAAsp), IAA-glutamate (IAGlu), 2-oxoindole-3-acetic acid (oxIAA), IAA-glucose and oxIAA-glucose) was monitored by mass spectrometry. Our results are in consistence with published role of particular enzyme in Arabidopsis metabolism[3]. Applied isolation procedure and LC-MS/MS method showed sufficient sensitivity for in vivo labelling studies. Therefore, our developed approach allows to monitor phytohormone metabolic pathways in extremely small samples (<1 mg), such as root tips, meristems or embryos.

Conclusions

Obtained results provide us with important information about the rates of IAA biosynthesis and degradation and the relevance of individual IAA precursors, catabolites and conjugates. Combined with mutant line screening[4], the developed approach can improve our knowledge about auxin metabolic pathways in different plant species.

Novel Aspect

Using heavy isotope labelling and mass spectrometry, we have developed a new functional bioanalytical

approach for studying IAA biosynthesis and metabolism in Arabidopsis.

References

1. Novák O., Hényková E., Sairanen I., Kowalczyk M., Pospíšil T., Ljung K., *Plant Journal*, 72, 523 (2012).
2. Novák O., Pěnčík A., Blahoušek O., Ljung K., *Current Protocols in Plant Biology* 1, 419 (2016).
3. Ljung K., *Development*, 140, 943 (2013).
4. Pěnčík A., Casanova-Sáez R., Pilařová V., Žukauskaitė A., Pinto R., Luis Micol J., Ljung K., Novák O., *Journal of Experimental Botany*, doi:10.1093/jxb/ery084 (2018).

Acknowledgment

This work was supported by the Croatian Science Foundation (project no. IP-2014-09-4359) and by the Ministry of Education, Youth and Sports of the Czech Republic (project no. NPU I LO1204).

TP-222 / CAN OBSTRUCTIVE SLEEP APNOEA BE DIAGNOSED BY REAL-TIME ANALYSIS OF EXHALED BREATH USING SECONDARY ELECTROSPRAY IONIZATION MASS SPECTROMETRY - A VALIDATION STUDY

Nora Nowak⁽¹⁾ - ***Anna Engler***⁽²⁾ - ***Damaris Schneeberger***⁽²⁾ - ***Sira Thiel***⁽²⁾ - ***Thomas Gaisl***⁽²⁾ - ***Martin Thomas Gaugg***⁽¹⁾ - ***Pablo Martinez-Lozano Sinues***⁽³⁾ - ***Anna-Sophie Stöberl***⁽²⁾ - ***Malcolm Kohler***⁽²⁾ - ***Renato Zenobi***⁽¹⁾

ETH Zurich, Department of Chemistry and Applied Biosciences, Zurich (1) - University Hospital Zurich, Department of Pulmonology, Zurich (2) - University Children's Hospital Basel, Department of Biomedical Engineering, Basel (3)

Keywords: breath analysis, secondary electrospray ionization, obstructive sleep apnea, biomarkers, metabolomics

Introduction

Since breath analysis shows a high potential as a non-invasive diagnostic tool, many explorative case-control studies investigating different diseases have been carried out in the recent years. However, only very few results have been validated in a larger cohort later on. Here, we present a blinded validation study for the diagnosis of obstructive sleep apnoea (OSA) by analyzing exhaled breath with secondary electrospray ionization mass spectrometry (SESI-MS).

Methods

Subjects with suspected OSA (n=150) were diagnosed by conventional in-laboratory respiratory polygraphy. In addition, they underwent exhaled breath analysis using SESI-MS. Metabolic breath patterns were analyzed and the diagnosis of OSA was predicted for the subjects using a panel of biomarkers found to be characteristic for the presence of OSA in our previous pilot study [1]. The person analyzing the data was blinded regarding clinical diagnosis.

Results

Between the pilot study and the validation study, some technical improvements on our SESI-MS setup were made. We upgraded to a commercial SESI source making the data challenging to compare. In a preliminary data analysis of 130 subjects measured so far, we were able to render the old and the new dataset comparable using empirical Bayes methods[2]. Thus, the data from the pilot study can be used for training the classification model for predicting the validation cohort.

Conclusions

A positive outcome of this study will strengthen real-time breath analysis by SESI-MS as a diagnostic tool tremendously and bring it a step closer to its application in clinical routine.

Novel Aspect

This is the first validation study for any biomarkers found in exhaled breath using secondary electrospray ionization mass spectrometry.

References

1. Schwarz EI et al., Thorax, 0, 1-8 (2015).
2. Johnson WE, Li C, Biostatistics, 8, 1, 118-127 (2007).

TP-223 / A DUAL SEPARATION/HIGH RESOLUTION ACCURATE MASS SPECTROMETRY SYSTEM FOR IMPROVED METABOLOME COVERAGE AND INCREASED THROUGHPUT

Ioanna Ntai⁽¹⁾ - *Martin Samonig*⁽²⁾ - *Aran Paulus*⁽¹⁾ - *Ralf Tautenhahn*⁽¹⁾ - *Amanda Souza*⁽¹⁾ - *Andreas Huhmer*⁽¹⁾

Thermo Fisher Scientific, Marketing, San Jose (1) - Thermo Fisher Scientific, Marketing, Germering (2)

Keywords: metabolomics, high-throughput, untargeted metabolomics, metabolite profiling, beer

Introduction

Complete coverage of a metabolome of interest can be challenging due to the diversity of small molecules. Here, we integrated reverse phase (RP) and hydrophilic interaction chromatography (HILIC) to maximize throughput and metabolome coverage. This dual separation system was coupled to an Orbitrap™ mass spectrometer and used for untargeted metabolomic analysis of beer samples, identifying compounds important in discriminating flavor and indicating spoilage.

Methods

Sixteen beers were purchased from retail stores. For each beer, samples were collected at time of purchase and after storage at room temperature or 4°C for 6 and 18 weeks. A Thermo Scientific™ Vanquish™ UHPLC system equipped with an autosampler with two injection units, was coupled to a Thermo Scientific™ Q Exactive™ HF mass spectrometer. Data were processed with Thermo Scientific™ Compound Discoverer software for unknown identification and differential analysis.

Results

We developed a dual liquid chromatography/mass spectrometry (LC/MS) system that combines RP and HILIC separation and we evaluated it for its reproducibility and separation capabilities of complex samples. Small polar metabolites from beer, like malate and succinate, were not retained during RP, but showed good reproducibility of retention times and MS signal during HILIC separation. Hydrophobic metabolites, such as apigenin and isoxanthohumol, were reproducibly retained only during RP. By combining the two complementary separation techniques, we were able to reproducibly analyze a larger number of metabolites.

Two independent UHPLC pumps were incorporated into the system to allow independent control of the two columns (RP and HILIC). This enabled the use of different solvents, additives and pH ranges and resulted in broader metabolite coverage. By equilibrating one column, while separation was carried out on the other, analysis time was decreased by 30% and the overall throughput of the method increased.

Conclusions

The increased throughput of the dual LC/MS system described above, did not compromise reproducibility or compound detection. It was successfully employed for the differential analysis of beer and provided robust indicators of flavor profiles and degradation. One can envision extending its application to any metabolomics application, where increased throughput is needed.

Novel Aspect

Innovative integration of RP and HILIC separations leads to improved metabolome coverage while simultaneously increasing throughput of untargeted metabolomic workflow.

TP-224 / INTELLIGENT MSⁿ WORKFLOW FOR IMPROVED METABOLOME COVERAGE AND INCREASED CONFIDENCE IN UNKNOWN IDENTIFICATION

Ioanna Ntai - Iman Mohtashemi - Ralf Tautenhahn - Graeme McAlister - Seema Sharma - Vlad Zabrouskov - Amanda Souza - Andreas Huhmer

Thermo Fisher Scientific, Marketing, San Jose

Keywords: metabolomics, compound identification, untargeted metabolomics, multi-stage fragmentation, data-dependent acquisition

Introduction

Compound identification is a bottleneck in untargeted metabolomics, hindering biological interpretation of results. Here, we describe a data-informed workflow that maximizes the number of metabolites interrogated by MS/MS and MSⁿ, while minimizing the acquisition of uninformative spectra. This workflow was used to analyze human plasma resulting in high confidence identifications, deeper metabolome coverage and enhanced biological knowledge generation.

Methods

Human plasma was purchased from NIST. Metabolites were extracted with methanol and injected on a Thermo Scientific™ Hypersil GOLD™ column. Instrumentation included a Thermo Scientific™ Vanquish™ UHPLC system and a Thermo Scientific™ Orbitrap Tribrid™ Mass Spectrometer with modified instrument control and data acquisition software. Data were analyzed using Thermo Scientific™ Mass Frontier software and Thermo Scientific™ Compound Discoverer™ software.

Results

During data-dependent MS/MS, ions are selected based on abundance, without any knowledge of biological relevance or type of ion. In a typical DDA experiment, we determined, that >40% of MS/MS spectra could be attributed to background ions. By enabling the automatic generation and implementation of a background exclusion list based on real-time feature detection in LC-MS data, background ion MS² spectra were practically eliminated (<0.1%), allowing for the analysis of more true sample components. Small molecules form different types of adducts and cluster ions during electrospray ionization. Highly abundant compounds may prevent the fragmentation of metabolites of lower abundance. By populating the inclusion list with the preferred ion for each metabolite, more compounds can be sampled by MS/MS and MSⁿ in a single run. Additionally, by automatically updating inter-run inclusion and exclusion lists during analysis, we can ensure that compounds not selected for MS/MS and MSⁿ will be prioritized during a subsequent injection.

Conclusions

The combination of MSⁿ and automatically generated inter-run inclusion and exclusion lists resulted in fragmentation of more unique metabolites and a greater number of metabolites confidently annotated. Application of this innovative workflow addresses the identification bottleneck of untargeted metabolomics studies and enables confident biological interpretation of the results.

Novel Aspect

Automated workflow for information-rich fragmentation data acquisition, designed to minimize irrelevant spectra and maximize metabolome coverage.

TP-225 / NON-TARGET METABOLOMIC ANALYSIS OF PEPERONI DI SENISE PEPPERS PGI (CAPSICUM ANNUUM L.) AND EVALUATION OF THEIR NUTRACEUTICAL PROPERTIES

Alberto Onzo⁽¹⁾ - **Donatella Coviello**⁽¹⁾ - **Giuliana Bianco**⁽¹⁾ - **Raffaella Pascale**⁽²⁾ - **Maria Assunta Acquavia**⁽¹⁾ - **Philippe Schmitt-Kopplin**⁽³⁾ - **Patrizia Iannece**⁽⁴⁾ - **Carmine Gaeta**⁽⁴⁾

Università degli Studi della Basilicata, Dipartimento di Scienze, Potenza (1) - Università degli Studi della Basilicata, Scuola di Ingegneria, Potenza (2) - Helmholtz Zentrum Munchen, Analytical BioGeoChemistry, Neuherberg (3) - Università degli Studi di Salerno, Dipartimento di Chimica, Fisciano (4)

Keywords: Capsicum Annuum L., Peperoni di Senise, FT-ICR-MS, Phytochemical Screening, hypoglycemic activity

Introduction

Peperoni di Senise peppers (*Capsicum Annuum* L.), a typical food product cultivated in Basilicata (Southern Italy), protected with a PGI quality mark, are known for their unique taste. A Mass Spectrometry-based phytochemical screening was performed to promote their functional food properties.

Methods

Direct injection negative-ion ESI-Ultra High Resolution Mass Spectrometry was used to obtain a metabolic fingerprint of the sample. High-resolution mass spectra were acquired on a Bruker (BrukerDaltonik GmbH, Bremen, Germany) solariX XR Fourier transform ion cyclotron resonance mass spectrometer (FT-ICR-MS) equipped with a 7T superconducting magnet and an ESI source. Moreover, antioxidant activity of the extracts was evaluated.

Results

UltraHigh Resolution Mass Spectrometry data obtained from analysis of the extract of Peperoni di Senise peppers were used to perform a rapid analysis of metabolome by converting accurate m/z values in putative elemental formulas in order to better understand the chemical composition of the sample. Molecular formula maps were obtained by making Van Krevelen diagrams, that lead to a direct visualization of different classes of metabolites. As expected, the most abundant compounds were polyphenols and carotenoids. A significant antioxidant activity was assessed also.

Conclusions

The results of the present study shed some light on metabolic composition of Peperoni di Senise peppers and highlighted their health-promoting properties.

Novel Aspect

This study improved the knowledge of metabolites occurring in Peperoni di Senise fruits and it can be used for comparative studies with other types of *Capsicum Annuum* L. fruits.

References

1. Wahyuni Y., Ballester A. R., Sudarmonowati E., Bino R. J., Bovy A. G., J. Nat. Prod. (2013) 76, 783–793.
2. Wahyuni Y., Ballester A. R., Tikunov Y., de Vos R. C. H., Pelgrom K. T. B., Maharijaya A., Sudarmonowati E., Bino R. J., Bovy A. G., Metabolomics (2013) 9, 130–144.
3. Di Sotto A., Vecchiato M., Abete L., Toniolo C., Giusti A. M., Mannina L., Locatelli M., Nicoletti M., Di Giacomo S., Journal of Functional Foods (2018) 40, 679-691.
4. Sricharoen P., Lamaiphan N., Patthawaro P., Limchoowong N., Techawongstien S., Chanthai S., Ultrasonics Sonochemistry 38 (2017), 629–639.

TP-226 / DETERMINATION OF METABOLIZATION OF NEW FLUORESCENTLY LABELED AUXIN-LIKE COMPOUNDS IN PLANT MATERIAL USING UHPLC-MS/MS

Barbora Pařízková⁽¹⁾ - *Asta Žukauskaitė*⁽²⁾ - *Thomas Vain*⁽³⁾ - *Peter Grones*⁽³⁾ - *Martin Kubeš*⁽⁴⁾ - *Martin Kieffer*⁽⁵⁾ - *Karel Doležal*⁽¹⁾ - *Stefan Kepinski*⁽⁵⁾ - *Miroslav Strnad*⁽²⁾ - *Stéphanie Robert*⁽⁶⁾ - *Ondřej Novák*⁽²⁾

Centre of the Region Haná for Biotechnological and Agricultural Research, Department of Chemical Biology and Genetics, Palacký University, Olomouc (1) - Centre of the Region Haná for Biotechnological and Agricultural Research, Laboratory of Growth Regulators, Palacký University and Institute of Experimental Botany, Academy of Sciences of the Czech Republic, Olomouc (2) - Umeå Plant Science Centre, Department of Forest Genetics and Plant Physiology, Swedish University of Agricultural Sciences, Umeå (3) - Centre of the Region Haná for Biotechnological and Agricultural Research, Department of Chemical Biology and Genetics, Palacký University, Olomouc (4) - Centre for Plant Sciences, University of Leeds, Leeds (5) - Umeå Plant Science Centre, Department of Forest Genetics and Plant Physiology, Swedish University of Agricultural Sciences, Olomouc (6)

Keywords: UPHLC/MSMS, auxin, fluorescence, in planta stability

Introduction

Auxin is an indole-derived plant hormone which significantly contributes to plant growth and development [1]. Using a fluorescent label to visualize auxin natural distribution in different physiological processes will provide new insights into auxin biology [2]. Here we present one-step purification and quantification methods used to determine the chemical stability of the new fluorescent derivatives and their potential in vivo metabolism in plant system.

Methods

The purification procedure based on liquid-liquid extraction was developed using hexane:water:methanol (1:1:1) as an extraction solution. After isolation step, an Acquity UPLC® I-Class System coupled to a triple quadrupole mass spectrometer Xevo™ TQ-S MS(Waters) was employed for detection of the analytes. Concentrations of all compounds were then calculated by an isotopic dilution method using the stable isotope labelled standards by deuterium and ¹³C.

Results

Applying the presented method, time-dependent uptake and metabolism dynamics of two new fluorescent auxin derivatives in roots of *Arabidopsis thaliana* seedlings was determined. After the treatment of 5-d old plants by 10 μM exogenous application of tested compounds, both fluorescent conjugates showed similar metabolism rate in vivo yielding 20-28 pmol/50 roots of free auxin in 3 h (5 % of their total content). Furthermore, the detected amount of free auxin analogue in short-term treatment (30 min) was close to the detection limit of the UHPLC-MS/MS method (in range 5.8-9.5 pmol/50 roots representing just 2% of the total content of fluorescent analogues). Despite a slight metabolism on the background, the evidences of biological activity of the whole fluorescent compounds were provided in various auxin-related bioassays.

Conclusions

A new detection method was developed to determine stability of two fluorescent analogues of plant hormone auxin in planta. Our analytical approach helped to reveal the time-dependent metabolism of new fluorescent auxin derivatives in plant tissues. Furthermore, the detected amount of free auxin analogue was proven not to be sufficient to cause the biological response.

Novel Aspect

A developed purification method is applicable for determination of a stability of various plant hormone derivatives and presents a useful analytical tool for plant research.

References

1. Enders T. A., Strader L. C., American Journal of Botany, 102, 180-196 (2016).
2. Hayashi K-i., et al., Proceedings of the National Academy of Sciences, 111, 1-6 (2014).

TP-227 / METABOLOMIC APPROACHES TO INVESTIGATE THE ROLE OF THE MITOCHONDRIAL REGULATOR ZC3H10 IN ADIPOCYTES

Silvia Pedretti - Matteo Audano - Emma De Fabiani - Donatella Caruso - Nico Mitro

Università degli Studi di Milano, Departement of Pharmacological and Biomolecular Sciences, Milan

Keywords: Zc3h10; LC-MS/MS; metabolism; mitochondria; adipocytes;

Introduction

Mitochondria play a crucial role in many cellular processes and they are essential organelles for the health of the cell. Beyond their contribution to energy production, they are key regulators of tissue development and cell differentiation. We recently isolated the new mitochondrial regulator zinc finger CCCH-type containing 10 (Zc3h10) [1] and here, we validated its role during adipocytes differentiation.

Methods

C3H/10T1/2 cell line (mesenchymal stem cells) can be differentiated into white adipocytes using a specific adipogenic cocktail. Quantification of different metabolites was performed with a liquid chromatography/tandem mass spectrometry (LC-MS/MS) on an API-4000 triple quadrupole mass spectrometer coupled with a HPLC system using a C18 column for amino acids and cyano-phase LUNA column for metabolites.

Results

Zc3h10 protein levels increases during C3H/10T1/2 differentiation. Zc3h10 silencing significantly affects adipocyte differentiation and mitochondrial activity. To demonstrate that Zc3h10 plays a role in energy metabolism we evaluated the intracellular levels of by-products belonging to the main metabolic pathways (*i.e.* glycolysis, tricarboxylic acid (TCA) cycle and pentose phosphate pathway (PPP), amino acids) using LC-MS/MS. Steady-state metabolomics indicated that lack of Zc3h10 led to decreased AMP, ADP, ATP and NADH levels. We also observed reduced levels of acetyl-CoA, α -ketoglutarate (α -KG), citrate, and succinyl-CoA and increased levels of glutamate. We also used metabolic tracers ([U-¹³C₆]-glucose, [U-¹³C₁₆]-palmitate or [U-¹³C₅]-glutamine) to confirm that the flow of energy substrates into the TCA is affected by Zc3h10 silencing.

Conclusions

Our results indicate that Zc3h10 expression increases during murine white adipocyte differentiation. Further, Zc3h10 silencing in white preadipocytes and adipocytes deeply impaired mitochondrial function, decreased adipogenic potential and altered metabolic profile.

Novel Aspect

These data annotate Zc3h10 as a new regulator of mitochondrial function and cell differentiation in adipocytes.

References

1. Audano M., Pedretti S., Cermenati G., Brioschi E., Diaferia G. R., Ghisletti S., Cuomo A., Bonaldi T., Salerno F., Mora M., Grigore L., Garlaschelli K., Baragetti A., Bonacina F., Catapano A. L., Norata G. D., Crestani M., Caruso D., Saez E., De Fabiani E., Mitro N., EMBO Reports, (2018) Mar 5, pii: e45531. doi: 10.15252/embr.201745531. [Epub ahead of print].

TP-228 / LC-MS/MS QUANTIFICATION OF THIOL-CONTAINING METABOLITES WITHIN THE DE NOVO GLUTATHIONE SYNTHESIS PATHWAY OVER THE LIFESPAN OF CAENORHABDITIS ELEGANS

Gavin Ferguson⁽¹⁾ - *Russell Pickford*⁽²⁾ - *Wallace Bridge*⁽¹⁾

University of New South Wales, School of Biotechnology and Biomolecular Sciences, Sydney (1) - University of New South Wales, Bioanalytical Mass Spectrometry Facility, Sydney (2)

Keywords: glutathione, thiol, quantification, LC-MS/MS, aging

Introduction

A consistent underlying index of aging is declining cellular levels of the tripeptide glutathione (GSH). This could be due to either limited cysteine availability or a decrease in the activity of the rate-limiting enzyme glutamate cysteine ligase (GCL). Insight into the availability of cysteine levels and the regulatory control of GCL will inform the development of therapeutic strategies aimed to overcome age-related dysregulation of the GSH synthesis pathway.

Methods

We have developed a rapid and sensitive LC-MS/MS method that can simultaneously quantify each thiol in the de novo GSH synthesis pathway (Cys, γ -GC and GSH together with GSSG) using a relatively unexplored thiol-derivatizing reagent, 4-(dimethylaminoethylaminosulfonyl)-7-chloro-2,1,3-benzoxadiazole (DAABD-Cl). This method was applied to a nematode model to determine if the dysregulation in GSH homeostasis is due to cysteine availability or GCL activity.

Results

The nematode, *Caenorhabditis elegans*, which has a median lifespan of approximately 17 days, presents a potentially useful model to investigate the age-related decline in GSH homeostasis, due to the presence of the metabolites and orthologues of the mammalian GSH synthesis pathway enzymes. Work to date has validated the LC-MS/MS method as having the dynamic range, sensitivity and selectivity to be suitable for the analysis of age-dependent changes of all the thiol and disulfide levels in nematode tissue samples. Using this method, thiol and disulfide levels were monitored in aged-matched wild-type cohorts. It was observed that γ -GC and GSH levels significantly declined over the first 10 days of adulthood. The GSH:GSSG molar ratio followed an age-related oxidizing shift, declining from 115:1 in day 1 old worms to 25:1 in day 10 old worms. Cysteine levels did not appear to show any age-dependent trends, indicating that factors other than the availability of cysteine are perhaps causing this age-related decline in GSH synthesis.

Conclusions

The age associated decline of the GSH pool is likely due to a progressive lowering of GCL activity. The corresponding decline in the GSH/GSSG ratio may be a function of a decline in the activity of the GSH recycling enzyme, glutathione reductase, or a depletion of its NADPH cofactor pool. Future work will expand the assay and quantify other key thiol metabolites to further our mechanistic understanding of the aging process in relation to glutathione homeostasis.

Novel Aspect

We developed a strategy using LC-MS/MS coupled with derivatization for the quantitative measurement of labile low molecular weight thiols in biological tissue.

TP-229 / INTEGRATION OF TANDEM MASS SPECTROMETRY MOLECULAR NETWORKING AND GAS-PHASE FRAGMENTATION REACTIONS FOR STRUCTURAL ANALYSIS OF FLAVONOID GLYCOCONJUGATES

Alan Cesar Pilon⁽¹⁾ - *Fausto Carnevale Neto*⁽¹⁾ - *Hawei Gu*⁽²⁾ - *Daniel Raftery*⁽²⁾ - *Vanderlan da Silva Bolzani*⁽³⁾ - *Norberto Peporine Lopes*⁽¹⁾ - *Ian Castro Gamboa*⁽³⁾

University of São Paulo, Faculty of Pharmaceutical Sciences, Ribeirão Preto (1) - University of Washington, Northwest Metabolomics Research Center, Seattle (2) - University of São Paulo State, Institute of Chemistry, Araraquara (3)

Keywords: Molecular Networking (GNPS), Dereplication, O-glycosyl-flavonoids, Gas-phase Fragmentation

Introduction

Gluconjugated flavonoids (GF) are important natural products in plants [1]. Their identification in the different isomeric forms, stereochemistry and positions still represent a great challenge. HPLC-MS/MS combined with bioinformatics proved to be efficient in analyzing and organizing large datasets of MS2. For this work, we integrate the Molecular Network (MS2) and the fragmentation patterns to the GF annotation in Chrysobalanaceae species.

Methods

Our approach uses massive molecular networks that embeds chemical - glycosylation type and position, nature of their aglycones, and the structure/linkage information of their glycan moieties - and biological information layers to highlight the flavonoid composition through structural-based organization. We exemplify this workflow by identifying the glycoconjugated flavonoids present in six Chrysobalanaceae species analyzed using LC-ESI-QToF-MS.

Results

The combination of Molecular Network capabilities, gas-phase fragmentation patterns with detailed taxonomical knowledge on flavonoid biochemistry allowed a rapid screening and annotation of several glyconjugate flavonoids including those with different glycosylation type and position, nature of their aglycones and the structure/linkage information of their glycan moieties. According to the literature, more than fifty flavonoids have been described, so far, from Chrysobalanaceae species, with a predominance of flavonol glycosides [2]. Based on the systematic inspection of all MN nodes related to the flavonoid groups, sixty-four metabolites were identified in the Chrysobalanaceae plant species, with a predominance of flavonol 3-O-glycosides. This is also the first report of methoxyl, O-di- and tri-O-glycosyl flavonols for this family.

Conclusions

The application of this strategy using LC-MS in automatic data-dependent acquisition (DDA) capability combined to molecular networking allowed the visualization and data comparison not only for known compounds but also provided information for the determination of correlated analogues of glycoconjugated flavonoids described for the first time in the Chrysobalanaceae family.

Novel Aspect

This is first time molecular networking, fragmentation pattern and chemotaxonomic information is combined to annotate glyconjugate flavonoids considering glycosylation type and position, aglycones and linkage information of their glycan moieties

References

1. Quideau S. *Angew. Chemis. Int.*, 45 (41), 6786-6787 (2006).
2. Carnevale Neto F., Pilon A. C., da Silva Bolzani V., Castro-Gamboa, I. *Phytochem. Rev.*, 12 (1), 121-146 (2013).

TP-230 / TYPICAL LC-MS METABOLOMICS WORKFLOW FOR PROFILING URINE SAMPLES OF PATIENTS WITH COLORECTAL CANCER

Ivan Plyushchenko⁽¹⁾ - *Timofey Bolotnik*⁽¹⁾ - *Dmitry Shakhmatov*⁽²⁾ - *Sergey Achkasov*⁽²⁾ - *Oleg Sushkov*⁽²⁾ - *Oleg Shpigun*⁽¹⁾ - *Igor Rodin*⁽¹⁾

Lomonosov Moscow State University, Chemistry Department, Moscow (1) - State Scientific Center of Coloproctology, Moscow (2)

Keywords: colorectal cancer, pattern recognition, metabolomics workflow, data processing methods

Introduction

As mentioned in recent research in USA about 50000 deaths have happened from colorectal cancer (CRC), also authors postulate the need to create a method for screening CRC in the early stages [1]. The comprehensive workflow for pattern recognition of three group CRC patients based on UPLC-MS with fully comprehensive algorithm for optimization methods of preprocessing and data acquisition were reported here.

Methods

All morning fasting urine samples (on an empty stomach, immediately after sleep) of control group volunteers (8 pcs.), colorectal cancer patients before surgical operation (20 pcs.) and after surgical operation (12 pcs.) were collected for this research.

The HPLC separation was conducted on a C18 column. The separation was carried out in a gradient elution mode with MS detection in a positive ion (TIC mode). "Dilute and shoot" technique was used for sample preparation.

Results

Three types of signal drift correction were used for reduction unwanted variations [2]: total ion current, quantile, median and without any method. ANOVA was applied, to test correction methods. Two methods for metabolites concentration normalization were examined [3]: mass spectrometry total useful signal, creatinine concentration and without normalization. Log transformation and three type of scaling methods (auto, mean, Pareto and without scaling) were tested.

So, the dataset was applied to assess all possible 48 permutations of normalization, correction, transformation and scaling followed by PLS-DA and sPLS-DA. Analyzing sPLS-DA result, we found that using a log transformation, regardless of the application of any other data processing methods, the first 10 features with the maximum values of their loadings were the same for all permutations. After an extraction from the entire set of data set only the top features, in all cases samples groups were fully resolved from each other by PCA, sPLS-DA, PLS-DA, RF, dendrograms and heat map.

Conclusions

We propose the pragmatic decision procedure for selection and evaluation of preprocessing methods in metabolomics studies. This algorithm is based on the estimating of the influence each stage using multivariate and univariate statistical analysis. Extraction selected features from raw data set will led to appropriate pattern recognition if the model is built right.

Novel Aspect

This LC-MS workflow with algorithm of evaluation, optimization and scheme for performing statistical analysis operations can be applied in other relatively short studies.

References

1. Siegel R. L., Miller K. D., Fedewa S. A., Ahnen D. J., Meester R. G. S., Barzi A., Jemal A. CA Cancer J Clin, 67, 177-193 (2017).
2. NezamiRanjbar M. R., Zhao Y., Tadesse M. G., Wang Y., Ransom H. W. Proteome Science, 11, 1-13 (2013).
3. Wu Y., Li L. Journal of Chromatography A, 1430, 80-95 (2016).

TP-231 / METABOLOMICS ANALYSIS OF URINE IN MICE WITH PROSTATE CANCER

Laura Covadonga Rodas Sánchez⁽¹⁾ - ***Pablo Rodriguez-González***^(1,2) - ***José Ignacio García Alonso***⁽¹⁾ - ***Pedro González-Menéndez***⁽²⁾ - ***Juan Carlos Mayo Barrallo***⁽²⁾ - ***Rosa M. Sainz Menéndez***⁽²⁾ - ***Ángeles López-González***⁽³⁾ - ***Coral Barbas***⁽³⁾

University of Oviedo, Department of Physical and Analytical Chemistry, Faculty of Chemistry, Oviedo (1) - University of Oviedo, University Institute of Oncology (IUOPA), Oviedo (2) - Universidad San Pablo CEU, Centre for Metabolomic and Bioanalysis (CEMBIO), Madrid (3)

Keywords: metabolomics, metabolite profiling, urine, mass spectrometry, cancer

Introduction

Cancer is one of the leading causes of morbidity and mortality worldwide and the number of new cases is expected to rise significantly in the future. Improvements in early detection and therapeutic evaluation are essential to decrease cancer mortality [1]. Metabolomics is increasingly applied for the identification of biomarkers for disease diagnosis, prognosis and risk prediction. The metabolic analysis in non-invasive biofluids such as urine may help to obtain a complete functional molecular picture of the biochemistry and ultimately the discovery of new cancer biomarkers [2].

Methods

The different properties of the huge diversity of chemical compounds mean that the complete metabolic analysis of urine is practically impossible with a single analytical method. Therefore, multiple analytical approaches have been used such as capillary electrophoresis-mass spectrometry (CE-MS), liquid chromatography-mass spectrometry (LC-MS) and gas chromatography-mass spectrometry GC-MS.

Results

A urinary metabolic fingerprinting analysis has been performed in an animal model in order to detect potential biomarkers that could help in the early detection of prostate cancer, which is the second most frequent type of cancer in men [3]. The comparison between healthy and diseased mice was carried out analyzing urine samples from Wild-Type (WT) mice and Transgenic Adenocarcinoma of the Mouse Prostate (TRAMP) mice of C57/BL6 strain. The experimental design offers a wide range of metabolites with a significant chemical diversity, due to the different platforms applied: very polar (CE-MS) non-volatile (LC-MS) or volatile metabolites (GC-MS). This strategy produces large amounts of complex datasets that must be handled with powerful computational analysis. An unsupervised method, such as principal component analysis (PCA), reduces the complexity of the data obtaining plots of easier interpretation. This approach is useful to see whether the two groups (TRAMP and WT) can be discriminated by their characteristics.

Conclusions

To summarize, we have a preliminary study that permits the detection of a wide range of metabolites, some of which could be possible biomarkers because it allows the discrimination of the tumour progression. Those compounds are identified based on the capabilities of mass spectrometry.

Novel Aspect

The search for biomarkers as an alternative to the analysis of the prostate-specific antigen (PSA) which has been questioned due to a lack of specificity at this cancer diagnosis [4].

References

1. O'Connell T, *Bioanalysis*, 4, 431-451 (2012)
2. Xia J., Broadhurst D., Wilson M., Wishart D., *Metabolomics*, 9, 280-299 (2013)
3. Struck-Lewicka W., Kordalewska M., Bujak, R., YumbaMpanga A., Markuszewski M., Jacyna J., Matuszewski M., Kaliszan R., Markuszewski M. J., *J. Pharm. Biomed. Anal.*, 111, 351-361 (2015)
4. Eastham J. A., Riedel E., Scardino P. T., Shike M., Fleisher M., Schatzkin A, Lanza E., Latkany L., Begg C. B., *JAMA*, 289, 2695-2700 (2003)

TP-232 / A NOVEL INTEGRATIVE STRATEGY TO PREVENT COLORECTAL CANCER WITHIN THE DIET-HOST-MICROBIOTA TRIANGLE: FROM ORGANIDS TO HUMAN IN VIVO REALITY

Josep Rubert⁽¹⁾ - ***Fulvio Mattivi***⁽²⁾ - ***Andrea Lunardi***⁽¹⁾

Centre for Integrative Biology (CIBIO), University of Trento, Trento, Italy (1) - Center Agriculture Food Environment (C3A), University of Trento, Trento, Italy (2)

Keywords: metabolomics, microbial metabolites, gut microbiota, cancer, organoids

Introduction

Currently, colorectal cancer (CRC) is the second cause of death in EU[1]. In EU, the Czech Republic and Slovakia showed the highest rates. By contrast, Greece presented the lowest rates across Europe[2, 3]. This incidence is postulated to reflect risk factors associated to socioeconomic status, including poor dietary habits and obesity. However, emerging evidence also implicates gut microbiota as an important effector in the relationship among diet, human health and cancer[4].

Methods

A novel integrative strategy to deeper understanding of CRC is vital. To accomplish this, a real-life scenario is needed. In this context, colon organoids/tumoroids will be established. Afterwards, an apple will be digested/fermented in vitro using a batch culture colonic model inoculated with feces from lean/obese healthy donors. In the end, such polyphenol metabolites will be tested in colon organoids/tumoroids and the mechanisms of action will be revealed by metabolomics and organoids assays.

Results

In summary, the results from TRIANGLE through the integration of polyphenols (diet), gut microbiota and colon organoids (host), will be highly multidisciplinary, and will undoubtedly set the stage towards the identification of mechanisms contributing to CRC understanding and will ultimately pave the way to phytochemical treatment and prevention. First, it is still unknown how apple phytochemicals are affected by lean and obese microbiotas, and it must be understood in order to comprehend chemopreventive efficacy of polyphenols and to be able to give nutritional recommendations related to functional-group class therein. Secondly, tumor organoids (tumoroids) and organoids can bridge the gap between human 2D cancer cell lines and animal-based models. Lastly, metabolomics analysis and 3D cell assays of colon organoids/tumoroids will provide valuable new insights into the mechanisms by which nutrient-gene interaction influences colon stem cell niche and CRC, and will open up new possibilities for CRC understanding and prevention.

Conclusions

TRIANGLE will provide novel insights towards CRC prevention by (1) using a new in vitro model based on colon organoids, (2) considering both gut microbiota composition and microbial metabolites, and (3) applying novel techniques. By implementing this, colon organoids will clearly mimic human real-life scenario. Secondly, the manner how polyphenols are affected by microbiotas will reveal the type and quantity of polyphenol metabolites, which could be correlated to the chemoprevention of CRC.

Novel Aspect

Metabolomics will provide a holistic signature in organoids and, together with organoid assay data can generate new hypotheses regarding underlying the mechanisms by which polyphenol-gene interaction influences CRC.

References

1. Torre L. A., Bray F., Siegel R. L. et al, CA Cancer J. Clin., 65, 87-108 (2015).
2. Ferlay J., Steliarova-Foucher E., Lortet-Tieulent J. et al, Eur. J. Cancer, 49, 1374-1403(2013).
3. Siegel R. L., Miller K. D., Jemal A., CA Cancer J. Clin., 66, 7-30 (2016).
4. Song M., Garrett W. S., Chan A. T., Gastroenterology, 148, 1244-1260 (2015).

TP-233 / INDOLOME ANALYSIS FOR NUTRACEUTICAL AND PHYSIOLOGICAL STUDIES

Federico Maria Rubino⁽¹⁾ - **Michele Dei Cas**⁽¹⁾ - **Jessica Rizzo**⁽¹⁾ - **Andrea Penoni**⁽²⁾ - **Marcello Iriti**⁽³⁾ - **Rita Clara Paroni**⁽¹⁾

Università degli Studi di Milano, Dipartimento di Scienze della Salute, Milano (1) - Università dell'Insubria, Dipartimento di Scienza e Alta Tecnologia, Como (2) - Università degli Studi di Milano, Dipartimento di Scienze Agrarie e Ambientali - Produzione, Territorio, Agroenergia, Milano (3)
 Indolome analysis for nutraceutical and physiological studies.

Keywords: isotope dilution, liquid chromatography, melatonin, tandem mass spectrometry, tryptophan ethyl ester

Introduction

Biotransformation products of tryptophan are gaining an increasing multi-faceted interest as nutritional supplements to cope with environmental and lifestyle related stressors. Their identification and measurement in nutritional and biological matrices mandates for flexible, yet accurate and precise analytical methods based on mass spectrometry and adapted to studies with stable isotopic labels. Fragmentation studies are necessary to fulfil this aim.

Methods

An API3000 LTQ LC-MS system is used for extensive tandem MS fragmentation studies on 3 tryptophan amino-acids, 4 tryptamines and 3 *N*-acetyl-tryptamines. A stable-isotope analogue of melatonin is custom synthesized. LC separation and quantification parameters are optimized for the detection by MRM of 8 analytes and 3 internal standards in four different biological and food matrices.

Results

Spectroscopic studies highlight key fragmentation pathways for the identification and measurement of 8 target analytes. The study of the fragmentation pattern [1] allowed designing and synthesizing [2] an under-considered isotope-labelled form of melatonin [3], to improve quantification by ID-MS-MS, with respect to previously employed isotopologues. Furthermore, the discrimination of the critical pair of isobaric melatonin and isomeric tryptophan ethyl ester, that also yield very similar fragmentation, could be accomplished with a scan function that addresses a triplet of consecutive fragments, only a pair of which are specific for each compound. With this method, the levels of melatonin and of its metabolic pathway are measured in human fluids from several unique conditions, as well as from nutraceutical and pharmaceutical formulations and food matrices.

Conclusions

Accurate fragmentation studies of biological indole compounds supply information to design a new internal standard for isotope dilution measurement of melatonin, and identification and quantification methods suitable for the measurement of the indolome in different natural matrices of interest for biological and technological applications.

Novel Aspect

An innovative internal standard to measure melatonin by ID-MS-MS and a method to discriminate isomers and isobars of melatonin are proposed.

References

1. Rubino F. M., Antognetti J., Campisi G., Rizzo J., Iriti M., Paroni R., Proc. 40th Meeting of the Mass Spectrometry Division of the Italian Chemical Society. Milano, Dec. 12 2014, 22-23 (2014).
2. Patent: Blagg J., Hercules Technology Management Co V, Inc. - WO2011/47156, 2011, A1
3. Hattox S. E., Murphy R. C., Biomed. Mass Spectrom. 2, 272 (1975).

TP-234 / RAPID ANALYSIS OF NCI60 PANEL METABOLIC AND LIPID PROFILES WITH AN AUTOMATIC WELL PLATE READER USING LASER ASSISTED REIMS

Richárd Schöffner⁽¹⁾ - *Julia Balog*⁽¹⁾ - *Daniel Simon*⁽¹⁾ - *Nora Kucsma*⁽²⁾ - *Anna Lovrics*⁽²⁾ - *Gergely Szakacs*⁽²⁾ - *Steven Pringle*⁽³⁾ - *Zoltán Takáts*⁽⁴⁾

Waters Research Center, R&D, Budapest (1) - MTA-TTK, Research Centre for Natural Sciences, Budapest (2)
- Waters Corporation, R&D, Wilmslow (3) - Imperial College London, South Kensington Campus, London (4)

Keywords: Cell line, metabolic profile, REIMS

Introduction

In this study, we present an automated, fast and effective way to identify different cell cultures based on their metabolic profile using laser assisted rapid evaporative ionization mass spectrometry (REIMS [1]). The aim of this study was to create an integrated workflow for the automatic analysis of cell line samples from well plates, and the building of a NCI60 cell line database based on the metabolic and lipid fingerprints of the samples.

Methods

The technology is based on the direct analysis of aerosol generated by a laser sampling device through thermal ablation of the samples. OPOTEK Opolette IR 532 tunable laser was combined with 3 stages (Thorlab NRT) for 3 dimensional sample movements. This stage prototype was coupled to a XEVO G2-XS QTOF MS and used on a number of different liquid and dried cell line samples from the NCI60 lines. Different multivariate statistical models were used to characterize cell lines.

Results

A total of 30 different cell lines were analyzed in both dried and liquid form in the mass-to-charge range of 50-2000. Number of parameters such as sample format (dried and liquid in different solvent, cell number per well), laser parameters (wavelength, pulse energy, frequency) and sampling methods (e.g. : aerosol intake) were optimized. The results showed that best signal was achieved with 2900 nm, 20Hz and a minimum of 3 mJ energy was applied on 1 Million cells / well. As the acquired mass spectra featured partly different species in liquid form compared to the dried samples, thus a different database was built for the two sample format after lock mass correction, background subtraction and normalization. Testing different classification models created with numerous artificial intelligence algorithms (Support Vector Machines, Neural Networks, Random Forests and Principal Component Analysis followed by Linear Discriminant Analysis) showed us a > 98.0% correct classification rate based on phospholipids, plasmalogens and sphingolipids.

Conclusions

The results demonstrate that our technology was capable of a rapid, fully automatic LA-REIMS based well plate reading and identification of different cancer cell lines. This method is suitable to analyze lipid and metabolic compounds of different cell lines. Our future plan is to finish the analysis of all cell lines within the NCI60 panel, and to identify the effect of perturbation of cell culturing environment on the metabolic and lipid profile of the cells.

Novel Aspect

Automated high throughput well plate analysis of NCI60 cell line metabolic fingerprint by OPO laser based ambient MS technique.

References

1. Strittmatter N., Lovrics A., Sessler J., McKenzie J. S., Bodai Zs., Doria M. L., Kucsma N., Szakacs G., Takats Z., *Analytical Chemistry*, 88, 7507-7514 (2016).

TP-235 / A WELLNESS STUDY USING MICROFLOW TARGETED METABOLOMICS TO INVESTIGATE THE EFFECTS OF DIET AND EXERCISE ON THE METABOLOME

Zuzana Demianova⁽¹⁾ - **Khatereh Motamedchaboki**⁽²⁾ - **Baljit K. Ubhi**⁽²⁾ - **Detlev Schleuder**⁽³⁾

SCIEX, Omics, Darmstadt (1) - SCIEX, Omics, Redwood City (2) - SCIEX Germany (3)

Keywords: metabolomics, QTRAP®, microLC, dried blood spot, polar metabolites

Introduction

LC-MS/MS is an essential tool for identification and quantitation of metabolites in complex samples to investigate the affected metabolic pathways in dietary assessment studies and precision medicine. Self-testing is on the rise as people are interested in monitoring their health. Here, we have analyzed dried blood spots (DBS) with microflow LC coupled to QTRAP® 6500+ system for a wellness study and monitored polar metabolites covering all major metabolic pathways.

Methods

Metabolite extracts from a 6-mm dry blood spot disk at three metabolic states; fasted, fed and active were analyzed in 5 replicates on a M3 MicroLC system with HILIC separation at pH9 using Luna- NH₂ columns (Phenomenex). MRM analysis performed in positive/negative polarity on a QTRAP® 6500+ system with IonDrive™ Turbo V source. 187 positive and 176 negative ion mode MRM's were combined into a single experiment to monitor polar metabolites across metabolic pathways.

Results

A total of 312 polar metabolites were monitored for this wellness study across the three metabolic states (fasted, fed and active). The PCA analysis of these groups in MarkerView™ Software clearly shows that these three groups are differentiated based on their metabolic profile. It is known that during exercise protein degradation occurs, releasing amino acids (AA) into the body including Tyr which is released from human skeletal muscle in the presence of insulin hence we observed the spike in Tyr concentrations in the samples after cycling activity was completed. This release of AAs from protein breakdown also releases branched chain AAs such as Ile and Leu, which are differentiating the active group. Technical replicate analyses of the DBS extracts shows there is insignificant variation between sampling which does not interfere with ability to differentiate these metabolic conditions. The method provides accuracy and precision across each of the different DBS extracts. DBS is a minimally invasive way of collecting blood samples with low cost.

Conclusions

The microflow targeted metabolomics analysis of DBS shows to be a rapid, sensitive, and accurate method for profiling polar metabolites in a variety of biofluids and could be utilized for monitoring these targeted polar metabolites in other wellness-based studies as the method successfully differentiated these three different conditions we monitored during this wellness study.

Novel Aspect

A sensitive microflow targeted metabolomics method allowing detection of over 300 polar metabolites across multiple biochemical pathways for studying wellness.

TP-236 / FECAL METABOLOMICS OF A MOUSE MODEL OF AUTISM

Emily Sekera - Troy Wood - Heather Rudolph - Stephen Carro

University at Buffalo, Department of Chemistry, Buffalo

Keywords: Metabolomics, MALDI, autism, biomarkers, FT-ICR

Introduction

Autism spectrum disorders (ASD) are a group of neurodevelopmental disorders lacking a clinical biomarker for diagnosis. In preliminary studies from our group, it was observed that the level of stercobilin, a metabolite produced in heme catabolism, is depleted in the fecal material of an animal model of autism. Here, we discuss creating an isotopologue standard for quantitation by MALDI FT-ICR MS.

Methods

Utilizing labelled stercobilin synthesized by adding 18-O water to our standard and heating it at 70°C for several hours, it is possible to create an isotopologue to quantify stercobilin levels. Fecal samples are chosen in pairs between the ASD model group and gender-matched controls. Stercobilin is then extracted and prepared for MS analysis. Analysis is completed by MALDI FT-ICR MS (n=5).

Results

Although the initial 18-O isotopologue of stercobilin used was shown to be stable through 20 days of storage, back labelling was observed for longer storage. By incubating the sample a second time under the same conditions, labelling efficiency of stercobilin was comparable but more labelling was observed at the labile sites with a lower amount of unlabeled material. Once a new batch of isotopologue was prepared, a new calibration curve was generated to quantify stercobilin within the fecal samples. It is from this work that an observed depletion of stercobilin was observed in the fecal material in a mouse model of ASD. An average depletion of 45% has been observed with a p-value of less than 0.001 (n=14). A less dramatic depletion was also observed with stercobilinogen; more samples will need to be analyzed to determine the statistical significance of the results.

Conclusions

The 45% depletion of stercobilin in the fecal material of an ASD mouse model relative to controls at a greater than 99% confidence level suggests that depletion of fecal stercobilin may serve as a potential ASD biomarker in humans. While less statistically significant, it appears that stercobilinogen, the precursor to stercobilin, is also depleted in the fecal samples of ASD mouse models.

Novel Aspect

Quantitative analysis of stercobilin, a putative autism biomarker, in mouse excrement utilizing isotopologue standards and MALDI FT-ICR MS.

TP-237 / LIQUID EXTRACTION SURFACE ANALYSIS AND DIRECT ESI/NANOESI MASS SPECTROMETRY FOR HIGH THROUGHPUT URINARY METABOLOMICS APPLIED TO MALARIA

Salah Abdelrazig - Catherine Ortori - David Barrett

University of Nottingham, School of Pharmacy, Nottingham, United Kingdom

Keywords: Liquid Extraction Surface Analysis, flow injection mass spectrometry, direct nanoESI mass spectrometry, urinary metabolomics, biomarker discovery

Introduction

Direct electrospray/nano-electrospray ionisation mass spectrometry (ESI/nanoESI-MS) have great potential as high-throughput approaches for untargeted urine metabolomics [1-5]. Here, we compare standard liquid chromatography-mass spectrometry (LC-MS) with ambient ionisation MS methods of liquid extraction surface analysis MS (LESA-MS), flow injection ESI-MS (FIE-MS) and nanoESI-chip-based infusion for clinical urinary metabolomics.

Methods

LESA-MS, chip-based infusion and FIE-MS using a high resolution-MS were developed and validated against a standard LC-MS method [6]. A mixture of 35 urinary metabolites were used to optimise and establish a high throughput sample preparation and data analysis protocols and the methods were then applied for biomarker discovery in malaria (n=88).

Results

The analysis of malaria/control(n=41/47) urine samples was validated using pooled QC samples interspaced

within the run. Malaria generated 7526 and 576 ions from LC-MS and FIE-MS, respectively. OPLS-DA showed clear separation and clustering of malaria from the controls with cross-validation of $R^2Y=0.810/0.993$, $Q^2=0.538/0.583$, sensitivity (80%) and specificity (77%) for LC-MS/FIE-MS, indicating comparable results. Altered levels of 30 and 17 metabolites were identified by LC-MS and FIE-MS, respectively, including pipercolic acid, taurine, 1,3-diacetylpropane, *N*-acetylspermidine and *N*-acetylputrescine and may have the potential of being used as biomarkers of malaria.

Conclusions

The developed direct ESI/nanoESI-MS demonstrated high throughput capability combined with the ability to differentiate between urine samples in disease and health state. LESA-MS requires no sample preparation and only involves placing urine samples onto a slide surface, gives the method a further credit of being very simple and easy to use. Therefore, these methods are recommended to be considered as a fast diagnostic/screening tool in clinical practice.

Novel Aspect

High-throughput surface analysis and direct ESI/nanoESI-MS methods for clinical biomarker discovery. Advanced data processing for direct ESI/nanoESI-MS.

References

1. Trunzer, M. ; Graf, D. ; Kiffe, M. Rapid Communications in Mass Spectrometry, 21, 937-944 (2007)
2. Froesch, M. ; Bindila, L. M. ; Baykut, G. ; Allen, M. ; Peter-Katalinic, J. ; Zamfir, A. D., Rapid Communications in Mass Spectrometry, 18(24), 3084-3092 (2004).
3. Zamfir, A. ; Vakhrushev, S. ; Sterling, A. ; Niebel, H. J. ; Allen, M. ; Peter-Katalinić, J., Analytical Chemistry, 76, 2046-2054 (2004).
4. Chetwynd, A. J. ; David, A. ; Hill, E. M. ; Abdul-Sada, A., Journal of Mass Spectrometry, 49, 1063-1069 (2014).
5. Chetwynd, A. J. ; Abdul-Sada, A. ; Hill, E. M., Analytical Chemistry, 87, 1158-1165 (2015).
6. Want, E. J. ; Wilson, I. D. ; Gika, H. ; Theodoridis, G. ; Plumb, R. S. ; Shockcor, J. ; Holmes, E. ; Nicholson, J. K., Nature Protocols, 5, 1005-1018 (2010).

TP-238 / URINARY METABOLIC PROFILE OF NEWBORNS WITH TRANSPOSITION OF GREAT ARTERIES UNDERGOING CARDIAC SURGERY WITH CARDIOPULMONARY BYPASS

Manuela Simonato⁽¹⁾ - ***Luca Vedovelli***⁽¹⁾ - ***Massimo Padalino***⁽²⁾ - ***Giovanni Stellin***⁽²⁾ - ***Stefania Carollo***⁽³⁾ - ***Andrea Bandini***⁽³⁾ - ***Virgilio Paolo Carnielli***⁽⁴⁾ - ***Paola Cogo***⁽⁵⁾

Istituto di Ricerca Pediatrica "Città della Speranza", PCare Laboratory, Padova (1) - Padova University Hospital, Pediatric and Congenital Cardiac Surgery Unit, Department of Cardiac, Thoracic and Vascular Sciences, Padova (2) - Padova University Hospital, Department of Medicine DIMED, Anesthesia and Resuscitation Institute, Padova (3) - Polytechnic University of Marche and Azienda Ospedaliero-Universitaria Ospedali Riuniti, Division of Neonatology, Department of Clinical Sciences, Ancona (4) - University of Udine, Department of Medicine, Udine (5)

Keywords: transposition of great arteries; cardiopulmonary bypass; untargeted metabolomics; newborns

Introduction

Complete transposition of the great arteries (TGA) is a common cyanotic congenital heart defect. Surgical correction, with the use of cardiopulmonary bypass (CPB), should be performed within the first month of life. CPB-assisted open heart surgery results in brain hypoperfusion and in a powerful systemic inflammatory response and oxidative stress [1,2].

Methods

Thirteen TGA newborns undergoing arterial switch operation with the use of CPB were recruited. Urine samples were collected before and after CPB and were analyzed using an untargeted metabolomic approach based on UHPLC-high resolution mass spectrometry.

Results

Since the post CPB metabolic spectra could be heavily “contaminated” by metabolites derived from administered drugs, we have constructed a library with drugs used during surgery and the relative metabolites retrieved from urine samples. This library was applied to our samples and a total of 2163 drug metabolites were correctly identified and excluded from the analysis. Afterward we detected over 14000 unique features in urine and 67 metabolites were significantly differentially abundant between pre and post CPB samples ($p < 0.02$; $FDR < 0.10$). Multivariate PCA classification model was established to differentiate between pre and post CPB urine samples.

Conclusions

This is the first report on the metabolic response to cardiac surgery in TGA newborns. Using UHPLC-high resolution mass spectrometry we identified metabolites whose concentrations appeared to be associated with clinical outcome. Metabolic profiling could be useful for patient stratification and directing investigations of clinical interventions.

Novel Aspect

We developed a high coverage metabolomic approach that has potential for biomarker discovery in CPB related damage.

References

1. Vedovelli L., Padalino M., D’Aronco S., Stellin G., Ori C., Carnielli V. P., Simonato M., Cogo P., Interactive Cardiovascular and Thoracic Surgery 24(3), 436 (2017).
2. Donmez A., Yurdakok O., Journal of Cardiothoracic and Vascular Anesthesia; 28(3), 778 (2014).

TP-239 / VOLATILOME ANALYSIS FOR PARKINSON'S DISEASE DIAGNOSTICS

Eleanor Sinclair⁽¹⁾ - *Drupad Trivedi*⁽¹⁾ - *Yun Xu*⁽¹⁾ - *Camilla Liscio*⁽²⁾ - *Phine Banks*⁽²⁾ - *Joy Milne*⁽¹⁾ - *Monty Silverdale*⁽³⁾ - *Tilo Kunath*⁽⁴⁾ - *Royston Goodacre*⁽¹⁾ - *Perdita Barran*⁽¹⁾

University of Manchester, Chemistry, Manchester (1) - Anatune, Analytical Chemistry, Cambridge (2) - University of Manchester, Neurology, Manchester (3) - University of Edinburgh, Biological Sciences, Edinburgh (4)

Keywords: Parkinson’s disease, biomarker, sebum, volatilome, super-smeller

Introduction

Parkinson’s disease (PD) is the 2nd most common age-related neurodegenerative disorder. The onset of symptoms such as tremor, slowed movements and muscle rigidity do not present in the early stages; however do form the basis of clinical diagnosis [1]. This often results in late diagnosis and treatment, highlighting an urgent need to identify a biomarker. Evidence from a super-smeller has shown there is a characteristic odour released from the skin of PD sufferers [2].

Methods

Sebum was collected non-invasively using gauze swabs from the upper back of 63 participants; 20 control subjects and 43 PD patients. Dynamic headspace was measured directly from swabs using thermal desorption-gas chromatography-mass spectrometry (TD-GC-MS). For three samples the instrument was

interfaced with an odour-port allowing our super-smeller to highlight regions of interest within the GC chromatogram.

Results

Headspace analysis was performed in two groups; firstly a discovery cohort consisting of 30 subjects was used to build a classification model in which metabolites of interest were identified. Partial least squares discriminant analysis was performed and validated with 5-fold cross validation alongside 26 permutation tests. The variables (n=17) which contributed to classification were selected using variable importance in projections (VIP) scores, where $VIP > 1$. The volatile components of the second group, a validation cohort of 30 different subjects, were targeted for the presence or absence of these identified biomarkers. Out of 17 metabolites, 9 were found within the validation dataset and subsequent statistical testing was performed. Further to this, regions of interest were identified by superimposing olfactogram data with the corresponding TD-GC-MS chromatograms. A 2 minute window associated with a strong PD-like scent in all 3 samples correlated to the elution time of 3 differentially expressed volatiles identified in both cohorts.

Conclusions

We demonstrate there is a characteristic volatile profile associated with the scent of Parkinson's disease. Several potential volatile biomarkers have been found to be statistically significant in the differentiation between disease and non-disease subjects and also associated to the PD-scent detected by a super-smeller.

Novel Aspect

A non-invasive method of sampling has enabled volatilome analysis of sebum produced from the skin of Parkinson's disease sufferers and produced a panel of putative biomarkers.

References

1. De Maagd G., Philip A., Pharm. Ther., 40, 504-532 (2015)
2. Morgan J., Lancet Neurol., 15, 138-139 (2016)

TP-240 / DETECTION OF METABOLIC SIGNATURES OF CROHN'S DISEASE USING AN INTELLIGENT MSⁿ-BASED DISCOVERY WORKFLOW

Ioanna Ntai - Amanda Souza - Ralf Tautenhahn - Gina Tan - Andreas Huhmer

Thermo Fisher Scientific, Chromatography and Mass Spectrometry, San Jose

Keywords: Crohn's disease, untargeted, metabolomics, biomarker, identification

Introduction

Crohn's disease is a chronic inflammation of the gastrointestinal tract where the etiology remains largely unknown. LC-MS profiling of serum metabolites has become a powerful tool for biomarker detection. Here, an untargeted metabolomics approach using high resolution accurate mass (HRAM) Orbitrap™ technology with dynamic acquisition for automatic data dependent MS/MS and multistage fragmentation was used to analyze serum from healthy and disease donors.

Methods

Metabolites were extracted using cold methanol from 18 serum samples consisting of healthy donors and disease donors with or without infliximab treatment. Supernatant was dried and re-suspended in acidified water and separated with a C18 column. Samples were analyzed by an Orbitrap Tribrid™ mass spectrometer with modified acquisition software. Differential analysis and compound annotation were applied using Thermo Scientific™ Compound Discoverer™ software.

Results

In this study, we explored the use of an intelligent LC-MSn untargeted metabolomics workflow for the identification of potential discriminant markers of Crohn's disease. We analyzed serum samples from three different populations.

Accurate m/z measurements determined molecular formula and putative annotation of thousands of metabolites by searching against ChemSpider. To increase annotation confidence, MS/MS data was searched against the mzCloud library. Multistage fragmentation (MSn) enabled structural elucidation of unknowns. The modified instrument control software automated inter-run inclusion and exclusion lists, facilitating MSn acquisition for more unique metabolites.

Differential analysis detected metabolite perturbations in the urea cycle and catabolism of amino acids in serum from disease donors when compared to that of healthy donors; these metabolic changes are often associated with inflammation. Some of those inflammation driven changes were minimized when disease donors were treated with infliximab.

Conclusions

The combination of an untargeted metabolomics profiling approach with intelligent high resolution MSn acquisition described here, addresses the identification bottleneck in untargeted metabolomics. This strategy presents a promising and facile workflow for the discovery and identification of novel disease biomarkers that could lead to further biochemical insights in disease progression and treatment outcome.

Novel Aspect

Combination of HRAM Orbitrap technology with automated MSn acquisition strategy for untargeted metabolomics enables detection of Crohn's disease metabolic signature.

TP-241 / MONOAMINE MAPPING BY MASS SPECTROMETRY IDENTIFIED BRAIN NUCLEI REGULATING ANXIETY IN A SEROTONIN DEFICIENCY MODEL

Yuki Sugiura - Eiji Sugiuyama

Keio University, School of Medicine, Tokyo

Keywords: Imaging MS, on tissue derivatization, neuroscience, monoamines

Introduction

In the brain, monoamines, dopamine (DA), norepinephrine (NE), and serotonin (5-hydroxytryptamine, 5-HT), act as neuromodulators, diffusing from the nuclei to affect slow-acting receptors of many neurons. Levels of monoamines are important not only for cognitive functions, including memory and learning, but also mood, appetite, and sleep. However, the relationship between actual monoamines' concentrations within the nuclei and behavioral output remains unclear. Here, utilizing imaging mass spectrometry (IMS), we generate a 3-dimensional mouse whole brain atlas of serotonin (5HT), dopamine (DA) and norepinephrine (NE), consisted of 30 coronal slice images.

Methods

After euthanasia, the mice brains were isolated, immediately frozen and sliced into 8 μm sections. 5HTd₄ solution was sprayed on them as an internal standard with an automated sprayer (SunCollect). To perform on-tissue derivatization of monoamines, DPP solution and then DHB solution were manually sprayed. Imaging were performed by a MALDI-FT-ICR MS (Bruker Daltonics). Three dimensional images was constructed using SCiLS lab software.

Results

Using the atlas, we unexpectedly found several brain nuclei co-contain 5HT and a catecholamine (DA or NE)

at high levels; among these, paraventricular nucleus of the thalamus (PVT) was 2nd 5HT rich nuclei in the brain. This novel observation suggests that 5HTergic system was integrated into catecholaminergic system by sharing several major nuclei.

Moreover, we attempted to reveal which 5HT-rich nuclei regulates mice behavior. Here, we used acute tryptophan deficient model (ATD) which exhibits decreased total brain 5HT followed by anxiety behavior. By generating and comparing two 3D monoamine-atlas of control and ATD-model brains, we successfully visualized specific reduction of 5-HT in PVT, which correlates anxiety-like behavior.

Furthermore, we show that the abnormal behavior and brain 5HT level were recovered by injecting [¹³C][¹⁵N]Trp into the ATD model. By visualizing [¹³C][¹⁵N]Trp conversion into [¹³C][¹⁵N]5HT, it is demonstrated that several brain nuclei including PVT have surprisingly high turnover of de novo 5-HT synthesis.

Conclusions

These results proved usefulness of brain-monoamine mapping by imaging MS to study relationship between monoamine concentrations within small brain nuclei and animal behavior.

Novel Aspect

3D brain atlas of monoamines created by imaging MS was utilized to study relationship between monoamine concentrations and animal behavior.

TP-242 / CONSTRUCTION AND APPLICATION OF A HIGH-RESOLUTION MS/MS LIBRARY INCLUDING RETENTION TIME INFORMATION FOR RAPID IDENTIFICATION OF ENDOGENOUS METABOLITES

Matthias Szesny⁽¹⁾ - *Shuang Zhao*⁽²⁾ - *Xian Luo*⁽²⁾ - *Wan Chan*⁽²⁾ - *Ulrike Schweiger Hufnagel*⁽¹⁾ - *Aiko Barsch*⁽¹⁾ - *Liang Li*⁽²⁾

Bruker Daltonik GmbH, Metabolomics, Bremen (1) - University of Alberta, Department of Chemistry, Edmonton (2)

Keywords: HMDB; metabolite identification; QTOF-MS; MS/MS library; RT correction

Introduction

High-resolution LC-MS is an important platform for metabolite detection and quantitation. However, for untargeted Metabolomics rapid, unambiguous and universal compound identification is still challenging. In this work, we report the construction of a library for relevant endogenous metabolite and its' successful application to human biofluids.

Methods

Standards were obtained from the Human Metabolome Database (HMDB). They were injected into an Intensity Solo C18 column via an Elute LC system and detected by a QTOF-MS (impact QTOF, all Bruker Daltonics) for acquiring MS/MS library spectra and the retention time determination (RT). For the analysis of biofluids the same setup was used following a standard operating protocol (SOP) for consistency.

Results

In this study a library from over 800 endogenous metabolites was created. It contains MS/MS spectra of 635 compounds acquired in positive mode and 474 negative mode spectra. Up to 5 collision energy levels were applied for each standard giving more than 6000 MS/MS spectra in total. For each metabolite, library fragment spectra were manually curated by confirming each fragment via a molecular formula. For unambiguous identification we determined the RT of each standard. Automatic metabolite identification was performed in the MetaboScape software based on matching of multiple parameters: precursor mass accuracy and isotopic pattern, RT, and MS/MS spectrum quality.

Conclusions

The established SOP and library was applied analyze to biofluids, *e.g.* plasma and urine, and metabolite identification results will be presented. Finally, we have examined the portability of this library for different instrumental conditions and for different labs.

Novel Aspect

Development of a high-confidence and rapid metabolite identification method by combining a high-resolution MS/MS library and RT information.

TP-243 / DEVELOPING A ROBUST LC-MS PANEL FOR ASSESSING RESIDUAL PLATELETS AND OTHER SAMPLE QUALITY ISSUES IN PLASMA

Warren K. Y. Tan^(1,2) - **Chester L. Drum**^(2,3,4)

NUS Graduate School for Integrative Sciences and Engineering, National University of Singapore, Singapore (1) - Cardiovascular Research Institute, National University Health System, Singapore (2) - Translational Laboratory in Genetic Medicine, Agency for Science, Technology and Research, Singapore (3) - Department of Surgery, Yong Loo Lin School of Medicine, National University of Singapore, Singapore (4)

Keywords: biomarker, pre-analytical variability, residual platelets, freeze-thaw, pre-centrifugation delay

Introduction

Plasma samples are commonly used for biomarker screening by LC-MS. However, pre-analytical variability in processing blood can lead to spurious results[1]. Our previous experience with thymosin beta-4 as biomarker interest us in the problem of residual platelets in plasma [2]. Proteomics-based markers however may not be conveniently incorporated into a metabolomics-only workflow, hence a panel of metabolic markers for sample quality control is ideal.

Methods

Alipidomic-focused screen for markers that is consistently higher in platelet rich compared to platelet poor plasma was conducted. Plasma was extracted using a single-phase method [3], then analyzed on a UPLC-QTOF system. Lipid identification and shortlisting was done using Lipid Data Analyzer 2.6.2 (LDA)[4]. A systematic literature search was also conducted to identify other reported pre-analytical sample quality issues and markers proposed for them.

Results

Molecular features 774.54 *m/z* and 746.51 *m/z* are the most consistently elevated in platelet rich plasma compared to platelet poor plasma. Base on tandem ms spectrum matching on LDA, these had been tentatively identified as plasmalogen phosphatidylethanolamine (P-PE) 36:4 and P-PE 38:4 respectively. Further validation are underway.

Common sample quality issues identified in the literature are: inappropriate use/non-use of anticoagulant, pre-centrifugation time delay, freeze-thaw cycles, hemolysis, sample handling, and sample storage. Markers for markers pre-centrifugation delays are the best studied: Sphingadienine-1-phosphate [5], ornithine/arginine ratio [6], ascorbic acid/lactic acid [7], cytokine levels [8], gamma-glutamyl-transferase and lactate dehydrogenase [9] have been proposed as markers. Gamma-glutamyl-transferase and lactate dehydrogenase [7] also marks for excessive freeze-thaw cycles. Hemolysis can be detected by hemoglobin or discoloration.

Conclusions

P-PE 36:4 and P-PE 38:4 are candidate markers for residual platelets in plasma samples. Because P-PE is a major constituent of cell membrane, we suspect it may be “housekeeping” in platelet membrane. We

intend to further characterize whether these markers can normalize for residual platelet content in plasma. Through validating the existing plasma quality markers in the literature, we hope to put together a robust LC-MS panel for assessing plasma sample quality.

Novel Aspect

Standardized protocols can reduce variability, however compliance may be infeasible. A LC-MS panel for plasma quality check would be transformative for biomarker screening.

References

1. V. V. Hernandez, C. Barbas, and D. Dudzik, "A review of blood sample handling and pre-processing for metabolomics studies," *Electrophoresis*, vol. 38, no. 18, pp. 2232–2241, 2017.
2. W. K. Y. Tan, K. Purnamawati, L. S. Pakkiri, S. H. Tan, X. Yang, M. Y. Chan, and C. L. Drum, "Sources of variability in quantifying circulating thymosin beta-4: literature review and recommendations.," *Expert Opin. Biol. Ther.*, vol.0, no.0, pp. 1–7, Mar. 2018.
3. Z. H. Alshehry, C. K. Barlow, J. M. Weir, Y. Zhou, M. J. McConville, and P. J. Meikle, "An efficient single phase method for the extraction of plasma lipids," *Metabolites*, vol. 5, no. 2, pp. 389–403, 2015.
4. J. Hartler, A. Triebel, A. Ziegl, M. Trötz Müller, G. N. Rechberger, O. A. Zeleznik, K. A. Zierler, F. Torta, A. Cazenave-Gassiot, M. R. Wenk, A. Fauland, C. E. Wheelock, A. M. Armando, O. Quehenberger, Q. Zhang, M. J. O. Wakelam, G. Haemmerle, F. Spener, H. C. Köfeler, and G. G. Thallinger, "Deciphering lipid structures based on platform-independent decision rules," *Nat. Methods*, vol. 14, no. 12, 2017.
5. X. Liu, M. Hoene, P. Yin, L. Fritsche, P. Plomgaard, J. S. Hansen, C. T. Nakas, A. M. Niess, J. Hudemann, M. Haap, M. Mendy, C. Weigert, X. Wang, A. Fritsche, A. Peter, H. -U. Häring, G. Xu, and R. Lehmann, "Quality Control of Serum and Plasma by Quantification of (4E,14Z)-Sphingadienine-C18-1-Phosphate Uncovers Common Preanalytical Errors During Handling of Whole Blood.," *Clin. Chem.*, vol.0, p. clinchem. 2017. 277905, Mar. 2018.
6. M. Jain, A. D. Kennedy, S. H. Elsea, and M. J. Miller, "Analytes related to erythrocyte metabolism are reliable biomarkers for preanalytical error due to delayed plasma processing in metabolomics studies," *Clin. Chim. Acta*, vol. 466, pp. 105–111, 2017.
7. J. P. Trezzi, A. Bulla, C. Bellora, M. Rose, P. Lescuyer, M. Kiehn topf, K. Hiller, and F. Betsou, "LacaScore: a novel plasma sample quality control tool based on ascorbic acid and lactic acid levels," *Metabolomics*, vol. 12, no. 6, pp. 1–11, 2016.
8. J. -E. Lee, J. -W. Kim, B. -G. Han, and S. -Y. Shin, "Impact of Whole-Blood Processing Conditions on Plasma and Serum Concentrations of Cytokines," *Biopreserv. Biobank.*, vol. 14, no. 1, p. bio. 2015.0059, 2016.
9. H. J. Kang, S. Y. Jeon, J. -S. Park, J. Y. Yun, H. N. Kil, W. K. Hong, M. -H. Lee, J. -W. Kim, J. Jeon, and B. G. Han, "Identification of Clinical Biomarkers for Pre-Analytical Quality Control of Blood Samples," *Biopreserv. Biobank.*, vol. 11, no. 2, pp. 94–100, 2013.

TP-244 / REVEALING THE HIDDEN PART OF BACTERIAL METABOLISM

Marion Thomas⁽¹⁾ - **Ekaterina Dairii**⁽²⁾ - **Pierre Loïc Saaidi**⁽²⁾ - **Alain Perret**⁽²⁾ - **Marcel Salanoubat**⁽²⁾

CEA - Genoscope, 2, Rue Gaston Crémieux, Evry (1) - CEA - Genoscope, Université Paris Saclay, Evry (2)

Keywords: Structural elucidation, metabolism, Mass Spectrometry, Nuclear Magnetic Resonance

Introduction

We studied the metabolic adaptation of the bacterium *Acinetobacter baylyi* ADP1 (ADP1) after a biotic stress that consisted in shifting the carbon source from succinate to quinate (1). Quinate metabolism in ADP1 has been investigated for decades (2). Yet, we spotted unexpected and unknown metabolites in quinate-grown cells(1). These orphan metabolites (OM) may participate in new pathways and represent entry points into the dark part of bacterial metabolism.

Methods

100 µg of a first metabolite of $m/z = 211$ (M211) was purified by preparative HPLC and its structure was analyzed by NMR. The interpretation of the ^1H , COZY, HSQC-HMBC (^1H - ^{13}C), ROESY, ^{13}C and ^1H - ^{15}N spectra, recorded on a 600 Hz NMR fitted with a cryoprobe.

Results

NMR analyses led us to conclude that the molecule was a (3-aminopropyl)-amine-hydroxybenzoic acid. Because of the small quantity of purified material, NMR data acquired did not allow us to unambiguously conclude on the position of the alcohol and carboxylic acid groups on the benzene ring of the molecule. To overcome this uncertainty, the compounds corresponding to the different six possible structures were chemically synthesized (Synthenova.com). Comparison of the NMR signals on one hand and high resolution (LTQ-Orbitrap Elite) LC/MS/MS data on the other hand between the purified metabolite and the reference standards led us to assign the structure of the compound of interest as 3 (3-aminopropyl) amine 4-hydroxybenzoic acid. This is a previously undescribed compound that is missing from all interviewed databases.

Conclusions

The identification of a second metabolite (M193) is under way. These structural elucidations are a first step in deciphering pathways that take place during quinate metabolism. We would next identify genes involved in their synthesis. We are using the complete library of deletion mutants of ADP1 (3) for determining their molecular phenotype (2600 metabolomes) and seek mutants in which OM are absent. In such mutants, the deleted genes thus participate in MO synthesis.

Novel Aspect

We solved the structure of a novel metabolite and investigate the genes involved in its biosynthesis for deciphering previously unreported metabolic pathways.

References

1. Stuani, L. et al. Novel metabolic features in ADP1 revealed by a multiomics approach. *Metabolomics* 10, 1223-1238 (2014).
2. Young DM, Parke D, Ornston LN. Opportunities for genetic investigation afforded by *Acinetobacter baylyi*, a nutritionally versatile bacterial species that is highly competent for natural transformation. *Annu. Rev. Microbiol.* 59, 519-51 (2005).
3. de Berardinis, V. et al. A complete collection of single-gene deletion mutants of *Acinetobacter baylyi* ADP1. *Molecular systems biology* 4, 174 (2008).

TP-245 / ASSESSMENT OF THE THS EXPOSITION IN ANIMAL MODELS BY MULTIPLATFORM METABOLOMICS ANALYSIS

Sonia Torres⁽¹⁾ - **Sara Samino**⁽¹⁾ - **Pere Rafols**⁽¹⁾ - **Neema Adhami**⁽²⁾ - **Manuela Martins-Greens**⁽²⁾ - **Xavier Correig**⁽³⁾ - **Noelia Ramirez**⁽⁴⁾

Rovira y Virgili University, Department of Electronic, Electrical and Automatic Control Engineering, Tarragona (1) - University of California, Department of Cell Biology and Neuroscience, Riverside (2) - Spanish Biomedical Research Centre in Diabetes and Associated Metabolic Disorders, Rovira y Virgili University, Barcelona (3) - IISPV, Rovira y Virgili University, Tarragona (4)

Keywords: LC/MS-based metabolomics, NMR-based metabolomics, MS imaging, multiplatform metabolomics, third handsmoke exposure

Introduction

Thirdhand tobacco smoke (THS) is a novel and poorly understood pathway of tobacco exposure produced

by the accumulation of secondhand smoke (SHS) on environmental surfaces that ages with time, becoming progressively more toxic[1]. Although the emerging evidence of THS harm, the molecular mechanisms of THS exposure are still unclear. Here, we present a multiplatform metabolomics approach for the characterization of metabolic disorders in liver of mice exposed to THS.

Methods

C57BL/6 mice were exposed for 24 weeks to THS, under conditions that mimic exposure of humans in homes of smokers. Aqueous and lipid liver extracts were analyzed by ^1H NMR and UHPLC-MS. AMIX and Chenomix programs were used to process NMR data and XCMS to LC-MS. Statistics and pathway enrichment were performed using the Metaboanalyst platform tools. Mass spectrometry images of liver were acquired by laser desorption-ionization MS, and processed with MSI software [2].

Results

Analysis of the lipid extracts by NMR confirmed that THS exposure increases the levels of triglycerides, polyunsaturated and monounsaturated fatty acids, which can lead to the development of fatty liver disease, as previously described [3]. UHPLC-MS analysis of the extracts allowed the identification of more than 30 triglycerides and 10 phosphocholines elevated in liver of THS-exposed mice. MS images showed lipid accumulation and increased lipid droplets due to THS-exposure. Statistically relevant ions in the images were further identified as 14 triglycerides. Untargeted metabolomics analysis of the aqueous extracts allowed the identification of 48 statistically significant metabolites. Pathways analysis showed that 19 of these metabolites are involved in the glutamate, glutathione and nicotine/nicotinamide metabolism pathways, which are related to nucleic acid and protein synthesis among others[4],[5],[6]. Furthermore, the dysregulation of these pathways confirms the oxidative stress nature of THS exposure health effects[7].

Conclusions

The combination of three different analytical platforms allowed the identification of more than 80 metabolites relevant for the health assessment of THS exposure. NMR provided reliable identifications and quantifications and UHPLC-MS was key for the identification of many dysregulated metabolites and the identification of lipids, whereas MSI allowed us to enrich the biological interpretation by the mapping of lipids in control and THS-exposed livers.

Novel Aspect

First metabolomics characterization of dysregulations in liver from THS-exposed mice, by MS and NMR multiplatform approach, including extracts analysis and tissue MS images.

References

1. M. Martins-Green et al., "Cigarette smoke toxins deposited on surfaces: Implications for human health," *PLoS One*, vol. 9, no. 1, pp. 1–12, 2014.
2. P. Ràfols et al., "RMSI: An R package for MS imaging data handling and visualization," *Bioinformatics*, vol. 33, no. 15, pp. 2427–2428, 2017.
3. N. H. Kim et al., "Association between cotinine-verified smoking status and risk of nonalcoholic fatty liver disease," *Liver Int.*, no. November 2017, pp. 1–8, 2018.
4. B. Hang et al., "Adverse health effects of thirdhand smoke: From cell to animal models," *Int. J. Mol. Sci.*, vol. 18, no. 5, 2017.
5. S. D. Yelamanchi et al., "A pathway map of glutamate metabolism," *J. Cell Commun. Signal.*, vol. 10, no. 1, pp. 69–75, 2016.
6. G. Wu, Y. Fang, S. Yang, J. R. Lupton, and N. D. Turner, "Recent Advances in Nutritional Sciences Glutathione Metabolism and Its Implications for Health 1," no. December 2003, pp. 489–492, 2018.
7. N. Adhami, S. R. Starck, C. Flores, and M. M. Green, "A health threat to bystanders living in the homes of smokers: How smoke toxins deposited on surfaces can cause insulin resistance," *PLoS One*, vol. 11, no. 3, pp. 1–19, 2016.

TP-246 / SIMULTANEOUS QUANTITATIVE ANALYSIS METHOD FOR HYDROXYPROLINE AND 4-HYDROXYGLUTAMIC ACID USING ISOTOPOLOGUE-SRM TO CLARIFY HYDROXYPROLINE METABOLISM

Makoto Tsuji

Daiichisankyo RD Novare Co., Ltd., Pharmaceutical and Biomedical Analysis, Tokyo

Keywords: 4-Hydroxyglutamic acid, HILIC, isotopologue-SRM, Hydroxyproline metabolism

Introduction

The main component of kidney stones is calcium oxalate. Oxalic acid is synthesized from hydroxyproline (HP) in vivo in addition to those from diet. We examined intermediate metabolites on HP metabolic pathway to clarify HP metabolism. 4-hydroxyglutamic acid (4HG) was found to be a marker of HP metabolism, and an simultaneous analytical method for HP and 4HG was established.

Methods

Culture media were acidified and deproteinized with acetonitrile and then centrifuged. The supernatant were injected into LC-MS/MS operated in ESI-positive mode. Chromatographic separation was achieved with Acquity BEH Amide column. HP and 4HG were quantified using $^2\text{H}_3$ -HP and $^2\text{H}_3$ -Lys as internal standards.

Results

Chromatographic separation of HP from other amino acids such as Met, Leu, Ile and Asp was achieved. The range of calibration for HP and 4HG were from 0.5 to 1000 microM and from 0.05 to 100 microM, respectively. Isotopologue-SRM (iSRM) [1] enabled to quantitate HP and 4HG simultaneously improving the linearity of calibration. 4HG concentration was elevated after adding HP to culture medium of human hepatocytes.

Conclusions

Although 4HG analysis method with derivatization have already been reported[2], the developed method is highly sensitive in addition to excellent operability not requiring derivatization. This method is considered useful for the study of HP metabolism.

Novel Aspect

Simultaneous quantitation of metabolite at lower concentration and substrate at higher concentration using iSRM.

References

1. Tsuji M. et al., J. Chromatogr. B, 953–954, 38 (2014).
2. Pitt JJ., Willis F., Tzanakos N., Belostotsky R., Frishberg Y., JIMD Reports, 15, 1 (2015).

TP-247 / AN UNTARGETED METABOLOMICS APPROACH TO GO BEYOND ASPIRIN CANONICAL EFFECT

Linda Turnu⁽¹⁾ - ***Chiara Maria Manega***⁽¹⁾ - ***Benedetta Porro***⁽¹⁾ - ***Alessandro Di Minno***⁽¹⁾ - ***Fabrizio Veglia***⁽²⁾ - ***Donatella Caruso***⁽³⁾ - ***Elena Tremoli***⁽¹⁾ - ***Viviana Cavalca***⁽¹⁾

Centro Cardiologico Monzino IRCCS, Unit of metabolomics and cellular biochemistry of atherotrombosis, Milano (1) - Centro Cardiologico Monzino IRCCS, Unit of biostatistic, Milano (2) - Università degli Studi di Milano, Dipartimento di Scienze Farmacologiche e Biomolecolari, Milano (3)

Keywords: metabolomics, untargeted analysis, aspirin, urine profile

Introduction

Pharmacometabolomics aims to highlight metabolic signatures after drug exposure. Although the canonical cyclooxygenase-1-mediated effect of Acetylsalicylic acid (ASA) is well defined, a substantial variability in drug response exists and the mechanisms at the basis of this variability are still poorly understood [1-2]. Applying an untargeted metabolomics approach, we aimed to reveal novel information about biochemical pathways modified by ASA treatment.

Methods

In this study we have developed an untargeted liquid chromatography-mass spectrometry method to define the urine metabolic profile of healthy subjects (n=7) before and 7 days after 100 mg once-daily ASA treatment.

Results

Through this untargeted method, we detected 2007 features: among them, 64 metabolites significantly differed ($p < 0.05$) after ASA assumption. Pathway analysis, performed on identified metabolites, revealed low levels of those metabolites involved in histidine, alanine, aspartate and glutamate and purine metabolisms, after ASA treatment. In addition, we observed the decrease of several short- and medium-chain acylcarnitines.

Conclusions

The data here reported reveal unusual pathways affected by ASA treatment, in particular, the decrease in urine acylcarnitines levels indicates an increase in fatty acid β -oxidation process. This observation may suggest non-canonical use of ASA in clinical situations characterized by energy depletion.

Novel Aspect

This study demonstrate the pivotal role of the untargeted metabolomics approach in exploring new clinical application of drugs.

References

1. Gum P. A., Kottle-Marchant K., Welch P. A., White J., Topol E. J. ; J Am CollCardiol, 41,961-5 (2009).
2. Hovens M. M., Snoep J. D., Eikenboom J. C., Van Der Bom J. G., Mertens B. J., Huisman M. V. ; Am Heart J, 153, 175-81 (2007).

TP-248 / EVALUATION OF HIGH SPEED, HIGH RESOLUTION DATA INDEPENDENT ACQUISITION FOR THE ANALYSIS OF METABOLOMIC FLUX, KINETICS AND PATHWAY MAPPING

Zuzana Demianova⁽¹⁾ - *Baljit K. Ubhi*⁽²⁾ - *Loren Olson*⁽²⁾ - *Jean-Baptiste Vincendet*⁽³⁾

SCIEX, Omics, Darmstadt (1) - SCIEX, Omics, Redwood City (2) - SCIEX France (3)

Keywords: metabolomics, flux, SWATH, data independent acquisition, quantitation

Introduction

Scientists across many disciplines are looking at metabolomics to answer challenging questions. Quantitative studies have focused on the steady-state metabolite levels where experimental variability can mask the true degree of metabolic regulation. Herein, we evaluate the merits of high resolution approaches to the measurement of the flux heavy isotopes in a metabolic pathway shows great potential to elucidate the regulation / kinetics of metabolic pathways.

Methods

MDCK cells were cultured under sterile conditions. A control and ¹³C labeled time course were created

(0,1,2,5,10 and 20 min). Controls was fed with normal media and flux samples were fed with stable isotope labeled $^{13}\text{C}_6$ -glucose media. Samples were washed 2x with 37C PBS lysed in icecold 70% MeOH. Cell supernatant was injected onto weak anion exchange chromatography. Acquisition was performed on TripleTOF 6600 in positive and negative mode with DDA or DIA.

Preliminary Data

Samples were analyzed by multiple means to create a library of observed metabolites. Careful creation of control samples (unlabeled) was important for the discovery of metabolites as well as the confirmation of species by mass accuracy, isotope ratio and MS/MS fragment pattern. Multivariate analysis (PCA and PLS-DA) of perturbed samples was particularly useful for identifying species that were undergoing significant flux of heavy atoms into a pathway. Major differences in the degree of kinetic flux were apparent between labeled controls and treatments that received pathway inhibitors.

imple unit resolution data dependent MS/MS, which is the most traditional style of LC/MS data collection, was employed for basic confirmation of metabolites. In addition, DIA and some targeted MS2 modes were also employed to evaluate their utility for kinetic measurements. The results support that the traditional DDA workflows were adequate for identifying and doing some basic quantitation of metabolites.

Conclusions

The stochastic nature of data dependent MS2 data collection limits the use of MS2 to identification only. Measurement of flux kinetics and identification of the locale of heavy atom incorporation is made possible by using targeted "scheduled" analyses as well DIA.

This work demonstrates the combining high levels of resolution and scan speed at the MS and MS2 level. Level of performance can generate information rich targeted and DIA datasets for metabolomics flux studies.

Novel Aspect

The use of data independent acquisition (vs. traditional DDA) to elucidate the structure and kinetics of metabolomics flux.

TP-249 / UTILIZING A STANDARDIZED AND TARGETED QUANTITATIVE ASSAY COUPLED TO HRAM MS FOR THE ANALYSIS OF SERUM METABOLITES OBTAINED FROM DIABETIC STUDY SUBJECTS

Anastasia Kalli⁽¹⁾ - *Sally Webb*⁽²⁾ - *Manuel Kratzke*⁽³⁾ - *Markus Langsdorf*⁽³⁾ - *Stephen Dearth*⁽³⁾ - *Hai Tuan Pham*⁽³⁾ - *Therese Koal*⁽³⁾ - *Andreas Huhmer*⁽²⁾

Thermo Fisher Scientific, CMD Life Science Research, Metabolomics, San Jose (1) - Thermo Fisher Scientific, CMD Life Sciences Research, San Jose (2) - Biocrates Life Sciences AG, Biocrates Life Sciences AG, Innsbruck (3)

Keywords: targeted metabolomics, standardization, quantitation, reproducibility

Introduction

Type 2 Diabetes, T2D, is a metabolic disorder characterized by decreased insulin sensitivity and abnormal hepatic glucose production. Metabolomics creates new capabilities for understanding metabolic disorders as it enables gaining insights into physiological and pathophysiological processes. We employed a standardized, quantitative method for targeted metabolic profiling of 408 metabolites in serum samples from T2D study subjects and healthy controls.

Methods

Sample preparation was performed with the AbsoluteIDQ® p400 HR Kit (Biocrates Life Sciences, AG). Briefly, 10 µl of plasma samples were pipetted onto the 96-well plate containing the internal standards.

Following derivatization, metabolites and internal standards were extracted and diluted for subsequent analysis. Mass spectrometric analysis was performed on a Vanquish™UHPLC system coupled to a Thermo Scientific™ Q Exactive™ HF Orbitrap mass spectrometer.

Results

We aim to monitor changes in metabolite concentrations in serum samples obtained from T2D study subjects and healthy controls. We have identified several metabolites that were significantly altered between the two study groups (p values $< 9.97E-04$). For instance, we observed increased concentrations of acylcarnitine AC(2:0), cholesterol ester CE(22:6), diglycerides DG(41:1, DG(36:4), sphingomyelins SM (36:1), SM(42:3), and phosphatidylcholines PC(38:6), PC(40:6) in the T2D serum samples versus the healthy controls. The concentration of serotonin, triclecyride TG(52:6), and phosphatidylcholines PC(40:4) were lower in the T2D serum samples. Acetylcarnitine AC(2:0), phosphatidylcholines PC(38:6), PC (36:6) and PC(40:6), serotonin and Ile, have been previously reported to be associated with higher likelihood of T2D, obesity, or they have shown differential changes in response to oral glucose tolerance test[1-4] in agreement with our findings. In terms of the method reproducibility, for 90% of all metabolites the % RSD was $< 20\%$.

Conclusions

By employing a standardized and quantitative analytical method and high resolution accurate mass (HRAM) mass spectrometry we were able to monitor metabolite changes in serum samples obtained from T2D study subjects and healthy controls. Overall, our results demonstrate that the method provides standardized and accurate quantitation for monitoring and quantifying changes between different study groups, while, providing high reproducibility

Novel Aspect

Combining a standardized analytical method with high resolution mass spectrometry for targeted metabolic profiling of T2D serum samples

References

1. Anna F., et. al., Diabetes, 62, 639-648, 2013
2. Adams S. H., et. al., J. Nutr. 139, 1073-1082, 2009
3. Yeon-Jung K., et. al., PLoS ONE, 11, 1-11, 2016
4. Ho E. H., et. al., Diabetes, 62, 2689-2698, 2013

Research use only. Not for use in diagnostic purpose

TP-250 / DEVELOPMENT OF A SIMULTANEOUS ANALYSIS OF THE FAT-SOLUBLE VITAMIN USING THE LITHIUM ION

Naohito Yamada - Takuya Matsui - Hideyuki Kuno

Japan Tobacco Inc., Toxicology Research Laboratories, Hadano

Keywords: Lithium ion, Fat-soluble vitamins, Simultaneous analysis

Introduction

In drug development, it is frequently demanded to investigate the fat-soluble vitamins when malabsorption of the fat-soluble vitamins is suspected as a side effect of treatment. Despite the need for a high sensitive analysis with a limited sample volume especially in a non-clinical study, a quantification of the fat-soluble vitamins using LC-ESI-MS/MS has been difficult because of the low ionization efficiency and fragmentation reaction. The aim of the work was to develop a high sensitive simultaneous analysis of the fat-soluble

vitamins in a biological sample using a post column infusion of a lithium ion.

Methods

Lithium ion was added post-column and the fat-soluble vitamins were analyzed by the multiple reactions monitoring (MRM) method. The transition of the MRM method used was m/z 286 (M)⁺ to 255, 300 (M)⁺ to 255, 384 (M)⁺ to 118, 400 (M)⁺ to 118, 423 ($M+Li$)⁺ to 387, 431 ($M+H$)⁺ to 165 and 451 ($M+H$)⁺ to 187 for vitamin A, retinoic acid, vitamin D3, 25OH-vitamin D3, 1,25(OH)₂-vitamin D3, vitamin E and vitamin K1, respectively.

Results

When vitamins A and D3 were ionized in the presence of the lithium ion, the (M)⁺ as well as ($M+Li$)⁺ was detected as a precursor ion. Furthermore, the MS/MS spectrums pattern obtained from (M)⁺ or ($M+Li$)⁺ significantly differed from the MS/MS spectrums pattern obtained from the ($M+H$)⁺. Surprisingly, when the standard sample of the fat-soluble vitamins was measured using a MRM method in the presence of the lithium ion, the intensity and S/N ratio of the peak of the vitamins A and D3 were much higher than those measured without lithium ion. The magnitude of their increase was 5 to 10 times.

Conclusions

In a simultaneous analysis of the fat-soluble vitamins, the ionization efficiency and fragmentation pattern were improved using a post column infusion of a lithium ion. The sensitivity and S/N ratio increased because the product ion from (M)⁺ or ($M+Li$)⁺ significantly differed from the product ion from ($M+H$)⁺.

Novel Aspect

A quantitative fat-soluble vitamins analysis with improved sensitivity and specificity in the presence of the lithium ion and detecting (M)⁺.

TP-251 / HUMAN METABOLITE IDENTIFICATION BY SEARCHING A HIGH QUALITY AND COMPREHENSIVE TANDEM MASS SPECTRAL LIBRARY

Xiaoyu Yang - Pedatsur Neta - Stephen E. Stein

Mass Spectrometry Data Center, NIST, Gaithersburg, MD, USA

Keywords: Tandem mass spectral library; Human metabolite identification; Data analysis; LC/MS/MS; In-source fragmentation.

Introduction

LC/MS/MS is a routine technique for compound identification in metabolomic studies. However, accurately identifying metabolites from millions of mass spectra is critical for data analysis. We applied an extended high quality and comprehensive reference tandem mass spectral library with >17,000 compounds, >141,000 precursors, and >671,000 spectra to identify human metabolites from LC/MS/MS data.

Methods

Mass spectra (MS₂, MS₃, MS₄) of each authentic compound were acquired in positive and negative modes at different collision energies for major ions including isotopic precursors and in-source fragment ions with Orbitrap Elite and Lumos instruments (HCD, FT-IT, IT) for extending a tandem mass spectral library. Human metabolites were identified by searching this library with MS Search program 2.3.

Results

The extended tandem mass spectral library [1, 2] contains 671,473 spectra, 141,927 precursor ions from 17,631 compounds. Of these, 4,005 are human metabolites. This library has been examined and applied in

metabolite identification in human plasma, urine, and milk samples that were run on an Orbitrap Fusion Lumos HCD mass spectrometer. For example, 325 metabolites and 154 in-source fragments from these metabolites by various neutral losses (*e.g.* $[M+H-NH_3]^+$ from tryptophan) were identified in a urine sample by searching the library. Over 700 human metabolites were identified with positive and negative precursor ions (*e.g.* $[M+H]^+$, $[M+Na]^+$, $[M+2H]^{2+}$, $[M-H]^-$). The in-source fragments were also identified in ~50% of the identified metabolites. Although the structures of the in-source fragments can be easily deduced from the original compounds, these sub-structures may not exist in the sample.

Conclusions

The high-quality reference tandem mass spectral library (17,631 compounds, 141,927 precursor ions, and 671,473 spectra) can be used for comprehensive metabolite identification. The spectra of in-source fragments can identify common fragmentation artifacts and high duplicate “unknown” spectra, confirm metabolite identification and help identify metabolites that are not in the library.

Novel Aspect

A high-quality reference tandem mass spectral library, including 4,005 human metabolites, facilitates accurate identification.

References

1. Yang, X., Neta, P., Stein, S., *J. Am. Soc. Mass Spectrom.*, 28 (11), 2280-2287 (2017).
2. Yang, X., Neta, P., Stein, S., *Anal. Chem.* 86 (13), 6393-6400 (2014).

TP-252 / TiO₂ ASSISTED LASER DESORPTION/IONIZATION MASS SPECTROMETRY FOR RAPID PROFILING OF CANDIDATE METABOLITE BIOMARKERS FROM ANTIMICROBIAL RESISTANT BACTERIA

Rutan Zhang⁽¹⁾ - **Qin Qin**⁽²⁾ - **Baohong Liu**⁽¹⁾ - **Liang Qiao**⁽¹⁾

Fudan University, Department of Chemistry, Shanghai, China (1) - Changhai Hospital, The Second Military Medical University, Changhai Hospital, Shanghai, China (2)

Keywords: TiO₂ P25, LDI-MS, Drug-resistant bacteria, Metabolomics, ROC analysis

Introduction

Antimicrobial resistance (AMR) is one of the most severe problems affecting public health and safety. (1) It is crucial to understand the mechanism of antimicrobial resistance from molecular level. Herein, we presented an innovative TiO₂ enhanced LDI method, which is able to recognize different types of bacteria at both species and strains level based on metabolites profiling by MVA analyses.

Methods

The workflow mainly involves sample preparation, data acquisition and MVA analyses. Briefly, P25 TiO₂ was utilized as LDI matrix. And then, a step by step procedure was applied with formic acid and 1-butanol acting as lysis solution together with salicylic acid and dopamine serving as charging additives in positive and negative mode respectively. Mass peaks were then aligned using an open source R package “MALDIquant” and normalized by total peak areas. (2)

Results

This TiO₂-LDI MS based metabolomics method can easily type four species of bacteria, including *E. coli*, *K. pneumoniae*, *P. aeruginosa* and *S. aureus*, based on PCA and PLS-DA models. The robustness for bacteria taxonomy was demonstrated by iterative validation using 48 species of clinical bacteria.

It is worth noting that this method also achieved fast bacterial AMR recognition, which possesses promising clinical diagnostic values. A multivariate ROC curve was obtained based on all obtained 11 potential biomarkers by Monte-Carlo cross validation (MCCV) (3) with an AUC of 0.955. Pseudo-unknown bacteria

test further confirmed AMR recognition capacity of our method.

In the end, pathway and metabolites network analysis were also performed to find responsible metabolites and mechanisms for AMR formation of ESBL+ *E. coli* at molecular level. Cysteine and methionine metabolism, purine metabolism, arginine and proline metabolism, glyoxylate and dicarboxylate metabolism, as well as peptidoglycan biosynthesis were found in disorder in ESBL+ *E. coli*.

Conclusions

This research presented an innovative TiO₂ LDI-MS method with good performance in the applications of bacteria identification, rapid AMR recognition and metabolomics analysis, which may provide an alternative to traditional LC-MS and/or NMR based metabolic analysis strategies.

Novel Aspect

This method presented a new metabolomics method based on LDI-TOF technique, which solved some crucial problems involving complexed background interferences, peaks alignment and normalization.

References

1. O'Neill, J., Review on Antimicrobial Resistance,(2016).
2. Gibb, S. ; Strimmer, K., *Bioinformatics*,28, 2270-2271 (2012).
3. Xia, J. ; Wishart, D. S., *Current protocols in bioinformatics*,55, 14.10.1-14.10.91 (2016).

TP-253 / METABOLIC CHARACTERIZATION OF THP-1 MACROPHAGES POLARIZATION USING LC-MS-BASED METABOLITE PROFILING

Alaa Abuawad⁽¹⁾ - ***David Barrett***⁽¹⁾ - ***Amir Ghaemmaghami***⁽²⁾ - ***Dong-Hyun Kim***⁽¹⁾

University of Nottingham, Advanced Materials and Healthcare Technologies, Nottingham, United Kingdom
(1) - University of Nottingham, Division of Immunology, Nottingham, United Kingdom (2)

Keywords: Macrophages, THP-1, metabolomics, LC-MS, polarization.

Introduction

Macrophages are heterogeneous and plastic immune cells that change their phenotype and function in response to environmental cues. M1 and M2 macrophages represent pro and anti-inflammatory phenotypes which show distinctively different functions [1, 2]. LC-MS metabolomics study was applied to characterize the metabolic profiles of M1 and M2 macrophages. This can be translated clinically to interrogate and manipulate macrophages phenotype.

Methods

THP-1 cells macrophages were polarized with either GC-MSF, LPS and IFN- γ or M-SCF and IL-4 to generate M1 or M2 polarised cells respectively. Unpolarised cells (M0) were set as a control. Then, LC-MS-based metabolite profiling was performed on an Accela system coupled to an Exactive(Thermo Fisher Scientific) operating with electrospray ionization (ESI), and ZIC-pHILIC(4.6 x 150mm and 5 μ m particle size, Merck Sequant) was used for chromatographic separation.

Results

The results showed that M1 and M2 macrophages have distinctive metabolic profiles. Clear separation was shown in OPLS-DA models between M0, M1 and M2 macrophages. Moreover, various metabolic pathways were significantly disturbed upon polarization including sphingolipids and pyrimidine metabolic pathways in M1 THP-1 macrophages. Whilst arginine and proline, alanine, aspartate and glutamate metabolisms, were the most significantly perturbed metabolic pathways in M2 THP-1 macrophages. Relevant key metabolites were successfully selected as potential biomarkers for each of functional phenotype, such as

[SP (22:0)] *N*-(docosanoyl)-sphing-4-enine, [SP (24:0)] *N*-(tetracosanoyl)-sphing-4-enine, [SP (16:0)] *N*-(hexadecanoyl)-sphing-4-enine, cytidine, cytosine, and 5-hydroxy tryptophan for M1 macrophages and *N*-acetyl glutamate, *N*-acetyl aspartate, and 4-Imidazolone-5-propanoate for M2 macrophages.

Conclusions

LC-MS metabolic profiling was used to characterize M1 and M2 macrophages. This approach efficiently demonstrated the distinct metabolic profile for M1 or M2 macrophages and revealed potential biomarkers for each phenotype. These biomarkers can be translated to explain the mechanisms of macrophages polarization and consequently, manipulate the cells phenotype. For example M2 is required to resolve autoimmune diseases, where M1 is active against tumor cells.

Novel Aspect

Discovery of biomarkers by LC-MS metabolomics can lead to understand polarization mechanism and then, development of targeted therapies for treating various diseases.

References

1. Genin M., Clement F., Fattaccioli A., Raes M., Michiels C., BMC Cancer, 15,577(2015).
2. Nicholson J. K., Lindon J. C., Holmes E., Xenobiotica, 29(11), 1181-9(1999).

TP-254 / INVESTIGATING THE DIFFERENCES IN THE MOLECULAR EFFECTS OF FREE AND NANOPARTICLE LOADED METHOTREXATE ON PROLIFERATIVE AND NONPROLIFERATIVE CANCER CELLS USING ADVANCE CELL BASED METABOLOMICS

Mohammad Al-Natour ⁽¹⁾ - Ali Alazzo ⁽¹⁾ - Amir Ghaemmaghami ⁽²⁾ - Dong-Hyun Kim ⁽³⁾ - Cameron Alexander ⁽¹⁾

University of Nottingham, School of Pharmacy, Division of Molecular Therapeutics and Formulation, Nottingham, United Kingdom (1) - University of Nottingham, School of Life Sciences, Division of Immunology, Nottingham, United Kingdom (2) - University of Nottingham, School of Pharmacy, Division of Advanced Materials and Healthcare Technologies, Nottingham, United Kingdom (3)

Keywords: Metabolomics, LC-MS, PLGA, Nanoparticles, Methotrexate

Introduction

Polymer nanomedicine has drawn the attention of researchers in recent years, especially drug loaded NPs¹. However, the impact of drug loaded NPs at the molecular and cellular levels could be different from the free drugs. These discrepancies might be attributed to different levels of uptake or cellular trafficking routes, yet it is not fully understood. Metabolomics is relevant tool for understanding such mechanisms involved in the response drug delivery systems.

Methods

Methotrexate (MTX) loaded and blank PLGA nanoparticles were fabricated using the conventional nanoprecipitation technique. The fabricated NPs in addition to free MTX were dosed to THP-1 macrophages and A549 lung epithelial cancer cells for 24 hours, after that the cellular metabolism was quenched and the intra cellular metabolites were extracted. Global LC-MS based metabolic profiling was used to study the metabolic changes resulted upon different treatments.

Results

THP-1 cells were sensitive to the NPs more than to MTX and most responses can be linked to the phagocytosis process, the NPs resulted in down regulation in TCA cycle and in glycolysis, at the same time redox homeostasis impairment was recorded. A549 cells were sensitive to MTX more than to nanoparticles and the responses metabolic responses indicated apoptosis. MTX has resulted in nucleotide metabolism

impairment, oxidative stress and shift in the energy metabolism towards fatty acid oxidation (FAO), this shift in energy metabolism leads to accumulation of acetyl CoA as end product of FAO, and because TCA cycle is malfunctioning, acetyl CoA is transformed to Acetoacetate which leads to ketoacidosis and apoptosis. The NPs did not have the same impact on A549 cells which might be due to the fact that the epithelial cells are not phagocytic. Being not proliferative might explain the resistance of activated THP-1 cells to MTX, especially because MTX is an antifolate agent that exerts its action during the S phase in proliferation process.

Conclusions

Global metabolic profiling approach was used to study the effects of MTX and MTX loaded PLGA NPs on THP-1 and A549 cells successfully. THP-1 cells were sensitive to the NPs more than to MTX and most responses can be linked to the phagocytosis process. A549 cells were sensitive to MTX more than to nanoparticles and the responses lead to apoptosis.

Novel Aspect

This work highlights the utility of metabolomics approach in understanding the difference in the molecular activity of a therapeutic agent when dosed as in free form or in NPs.

References

1. Vicent MJ, Ringsdorf H, Duncan R., *Advanced Drug Delivery Reviews*, 61, 1117-1120(2009)
2. Govender, Thirumala, Stolnik, Snjezana, Garnett, Martin C., Illum, Lisbeth, Davis, Stanley S., *Journal of Controlled Release*, 57, 171-185 (1999)

TP-255 / CHARACTERIZATION OF PHENOLIC COMPOUNDS AND DETERMINATION OF ANTI-INFLAMMATORY AND ANTI-CHOLINESTERASE ACTIVITY IN 2 COLOMBIAN FRUITS: BOROJOA PATINOI AND BACTRIS GASIPAES

Daniel Arias Ramirez

Universidad de los Andes, Chemistry Department, Bogota, Colombia

Keywords: Borojoa patinoi, Bactris gasipaes, HPLC-PDA-ESI-MSn, polyphenols

Introduction

The endemic Colombian fruits are an abundant source of secondary metabolites especially of phenolic compounds [1], [2]. Borojoa patinoi and Bactris gasipaes, two fruits from the region, are possible sources of this type of metabolites. Besides, there isn't any research focused on the identification of phenolics in these two fruits. The aim of this research is the total quantification and identification of phenolic compounds in these two fruits.

Methods

The samples were acquired at Buenaventura's City market located at the Pacific coast of Colombia. The extraction was performed using methanol acidified with formic acid [3], [4]. The separation was carried out using a Kinetex C18 column (150 mm x 3 mm i.d., 3 μ m) [3]. The mobile phase used was 0.1% formic acid in water and methanol. The Ion Trap was operated in data-dependent, full scan and MSn mode to obtain fragment ion m/z with a collision energy of 35% and an isolation width of 2 m/z .

Results

The extraction of polyphenols was previously optimized for HPLC-ESI-MS/MS analysis. The extraction was carried out according to the obtained profiles from the methanolic extracts allowing the identification of new compounds by means of data bases [5] [6]. In addition, by means of the Folin-Ciocalteu method the total content of phenolic compounds was quantified in these extracts. Additionally, the antioxidant activity

was measured by the ABTS and DPPH methods. The anti-inflammatory activity by FOX method [7] and the anti-cholinesterasic activity by de Ellman method [8].

Conclusions

At the structures elucidation ten new compounds were found for *Bactris gasipaes* among which 8 polyphenols stand out. In *Borojoa patinoi* new compounds were found eight of which are polyphenols. IC₅₀/EC₅₀ values were also found using quercetin as a positive control for antioxidant activities and anti-inflammatory assay. EC₅₀ values were determined for the anticholinesterasic activity using Eserine as positive control.

Novel Aspect

This is the first study focused on the search for polyphenols for these two Colombian fruits

References

1. Osorio, C. & Almanza, O. Antioxidant activity of anthocyanin-rich colombian tropical fruits. in ACS Symposium Series (2013).
2. Contreras-Calderón, J., Calderón-Jaimes, L., Guerra-Hernández, E. & García-Villanova, B. Antioxidant capacity, phenolic content and vitamin C in pulp, peel and seed from 24 exotic fruits from Colombia. Food Res. Int. (2011).
3. Mascherpa, D., Carazzone, C., Marrubini, G., Gazzani, G. & Papetti, A. Identification of phenolic constituents in cichorium endivia Var. crispum and Var. latifolium salads by high-performance liquid chromatography with diode array detection and electrospray ionization tandem mass spectrometry. J. Agric. Food Chem. 60, 12142–12150 (2012).
4. Grotewold, E. The science of flavonoids. The Science of Flavonoids (2006). 28822-2
5. Food Database. (n. d.). FoodDB. [online] Available at: <http://foodb.ca/> [Accessed 23 May 2018].
6. Allen, F., Greiner, R. and Wishart, D. (2018). Competitive fragmentation modeling of ESI-MS/MS spectra for putative metabolite identification. [online] CFM-ID. Available at: <http://cfmid.wishartlab.com/> [Accessed 31 May 2018].
7. Pinto, M. C., Tejada, A., Duque, A. L. & Macias, P. Determination of lipoxygenase activity in plant extracts using a modified ferrous oxidation. Xylenol orange assay. J. Agric. Food Chem. 55, 5956–5959 (2007).
8. Nicola, C., Salvador, M., Escalona Gower, A., Moura, S. & Echeverrigaray, S. Chemical constituents antioxidant and anticholinesterasic activity of *tabernaemontana catharinensis*. Sci. World J. (2013).

TP-256 / PHENOLICS, ALKALOIDS AND LIPIDS OF *OCHROMA PYRAMIDALE*, *TAPIRIRA GUIANENSIS* AND *INGA EDULIS*, HOST PLANTS OF THE CATERPILLAR *LONOMIA DESCIMONI*

Felipe Benavides - Chiara Carazzone

Universidad de los Andes, Chemistry Department, Bogotá D. C., Colombia

Keywords: *Lonomia*, Host-plant, Secondary Metabolites, Primary Metabolites

Introduction

The *Lonomia* caterpillars are venomous animals found in Latin America and Colombia owns a great variety species of this genus. [1][2] Very few works has been realized about the *Lonomia* caterpillars in Colombia and their host plants. The aim of this investigation is to study of the polyphenols, alkaloids, lipids and proteins content of *Ochroma pyramidale*, *Tapirira guianensis* and *Inga edulis*, host plants of the caterpillar *Lonomia descimoni*.

Methods

The specimens were collected in the Colombian Amazon jungle. The leaves were macerated with liquid nitrogen and the extraction was performed with different solvents depending on the metabolites of

interest (polyphenols[3], alkaloids[4], lipids[5] and proteins[6]). The extracts were filtered and injected onto HPLC-ESI-MS/MS system. The spectra and chromatograms were analyzed for the identification of the metabolites present in the extracts.

Results

The extraction conditions for the lipids[3] and the proteins[4] were optimized previously to the analysis by HPLC-ESI-MS/MS. The alkaloids and polyphenols extraction conditions were performed according to [5] and [6], respectively. The profiles obtained of the different extracts of the leaves of the three plants under study allowed the comparison of the metabolites present. The tandem mass spectra of the metabolites were assigned and, when possible, identified. Additionally, we could establish that despite the presence of some metabolites with the same m/z relation in the three plants, these corresponded to different compounds, since the MS/MS fragments of these were different.

Conclusions

The chromatograms and the MS / MS spectra allowed us to observe differences between the polyphenol, alkaloid, lipid and protein profiles of the three plants under study. Only a few compounds are present in all plant species, therefore it would be interesting to investigate the effects of these metabolites on the diet of the caterpillar.

Novel Aspect

In spite of the Colombian biodiversity, most species remain uninvestigated. This is the first study on the host plants of *Lonomia descimoni*, a venomous caterpillar.

References

1. Carrijo-carvalho, L., & Chudzinski-tavassi, A, *Toxicon*, 49, 741–757, 2007.
2. Pineda, D., Amarillo, A., Becerra, J., Montenegro, G. *Biomédica*, 21, 328–332, 2001.
3. Carazzone, C., Mascherpa, D., Gazzani, G., & Papetti, A. *Food Chemistry*, 138, 1062–1071, 2013.
4. Yang, Y. -F., Qu, S. -J., Xiao, K., Jiang, S. -H., Tan, J. -J., Tan, C. -H., & Zhu, D. -Y. *Journal of Asian Natural Products Research*, 12, 1005–1009, 2010.
5. Silitonga, A., Masjuki, H., Ong, H., Mahlia, T., Kusumo, F, *Energy Sources, Part A: Recovery, Utilization and Environmental Effects*, 39, 1167-1175, 2017.
6. Daim, L., Ooi, T., Yusof, H., Majid, N., Karsani, S. *Protein Journal*, 34, 304–312, 2015.

TP-257 / MULTI-OMICS PROFILING OF PANCREATIC CANCER STEM-LIKE CELLS

Claudia Di Carlo⁽¹⁾ - ***Marcello Manfredi***^(2,3) - ***Jessica Brandi***⁽¹⁾ - ***Ilaria Dando***⁽⁴⁾ - ***Elisa Dalla Pozza***⁽⁴⁾ - ***Arianna Buzzi***⁽³⁾ - ***Emilio Marengo***⁽³⁾ - ***Marta Palmieri***⁽⁴⁾ - ***Daniela Cecconi***⁽¹⁾

Department of Biotechnology, Proteomics and Mass Spectrometry Laboratory, University of Verona, Verona, Italy (1) - ISALIT, Spin-off of Department of Sciences and Technological Innovation, University of Piemonte Orientale, Alessandria, Italy (2) - Department of Sciences and Technological Innovation, University of Piemonte Orientale, Alessandria, Italy (3) -Department of Neuroscience, Biomedicine and Movement, University of Verona, Verona, Italy (4)

Keywords: pancreatic cancer stem cells; pancreatic cancer; metabolomic and proteomic approaches; pancreas; deregulated pathways.

Introduction

Pancreatic ductal adenocarcinoma (PDAC), one of the most aggressive solid tumours of the pancreas, is characterized by a remarkable resistance that also depends on the presence of the so-called cancer stem cells (CSCs) [1,2]. It is evident that, to make an impact on pancreatic cancer, it is necessary to eradicate

CSCs with targeted therapies and, for this reason, a complete molecular characterization of CSC biology is fundamental.

Methods

The MIA PaCa-2, PaCa3 and Panc-1 CSCs and relatives parental cells were cultivated and subjected to metabolomics and proteomic analyses. Different protocols have been evaluated to improve the identification of metabolites from adherent and suspended cell lines. Different extraction approaches and different derivatization approaches have been evaluated and compared.

Results

The results obtained indicate that CSCs are characterized by downregulation of succinate, fumarate, and malate and by upregulation of glutamate, in comparison to the parental cell lines. A reduction of Krebs cycle intermediates and an increase of a Krebs cycle-related metabolite may suggest a switch to a less oxidative metabolism. The results of the proteomic analyses are still in the collection phase. Preliminary data suggest that a total of 121, 188 and 186 proteins are modulated in Panc-1, PaCa-3 and MIA PaCa-2 CSCs as compared to the parental cells with p-value ≤ 0.05 and fold change ≥ 1.5 , respectively.

Conclusions

Altogether our findings may clarify some critical aspects of the metabolic network signature of PDAC and may serve as a relevant pool of metabolites and proteins from which biomarkers and therapeutic targets can be identified. Further analysis of biological replicates and validation of results are required to confirm our preliminary data.

References

1. Dalla Pozza E., Dando I., Biondani G., Brandi J., Costanzo C., Zoratti E., Fassan M., Boschi F., Melisi D., Cecconi D., Scupoli M. T., Scarpa A. & Palmieri M. *International Journal of Oncology*, 46, 1099-108 (2015).
2. Brandi J., Dando I., Dalla Pozza E., Biondani G., Jenkins R., Elliott V., Park K., Fanelli G., Zolla L., Costello E., Scarpa A., Cecconi D. & Palmieri M. *Journal of Proteome Research*, 15, 310–322 (2017).

TP-258 / INFORMED UNTARGETED METABOLOMICS USING LC-MS/MS

Julica Folberth⁽¹⁾ - ***Olaf Jöhren***⁽²⁾ - ***Markus Schwaninger***⁽¹⁾ - ***Alaa Othman***⁽²⁾

University of Luebeck, Institute for Experimental and Clinical Pharmacology and Toxicology, Luebeck, Germany (1) - University of Luebeck, Bioanalytic Core Facility, Center of Brain, Behavior and Metabolism (CBBM), Luebeck, Germany (2)

Keywords: untargeted metabolomics, data-dependent acquisition, LC-MS/MS, Orbitrap, library

Introduction

Metabolite identification is the main challenge in untargeted LC-MS based metabolomics. Although there are libraries available, they suffer some limitation. Therefore an LC-MS² standards library of high resolution and accurate mass spectra was generated. Another aspect is the inevitable bias in every metabolomics study, due to the diversity of the metabolome. To assess this bias the LC-MS² standards library can be used (Informed unTargeted Metabolomics).

Methods

An automated flow injection/t-SIM/ddMS² method was developed to generate various MS² spectra of more than 600 synthetic standards of metabolites (MS² library). Chromatographic data was acquired under different conditions (LC library). Using a hybrid quadrupole-Orbitrap mass spectrometer (Q-Exactive™, Thermo scientific™) allows high-resolution and accurate-mass detection. Data analysis was performed by

using FreeStyle, Library Manager and Tracefinder.

Results

We developed a method to acquire MS² spectra at different collision energies without the need for manual injection. Data-dependent acquisition proved superior to other scan modes as it offers high-quality MS¹ and MS² spectra. Due to flow injection, several different MS² spectra could be acquired for every metabolite. Using this flow injection/t-SIM/ddMS² method, we generated about 4000 MS² spectra for over 500 metabolites. LC data was acquired for 295 metabolites in 139 unique conditions. Chemical properties of the metabolites will be included in the LC-MS² library as well. To further expand the library, more compounds, LC conditions and also extraction methods will be tested. By utilizing a certain extraction and LC method, an analysis is shifted to a subset of metabolites, that can be estimated by comparing features of metabolites already detected using this workflow and documented in the library. Informed unTargeted Metabolomics can be adjusted to cover the highest number of metabolites possible or to analyze specific metabolites of interest.

Conclusions

For unambiguous metabolite identification high mass accuracy, high-quality MS² spectra and retention times under specific chromatographic conditions are essential – information that is included in the LC-MS² standards library together with chemical properties of the metabolites. Informed unTargeted Metabolomics uses this information to assess or determine the bias that goes along with the impossibility to cover the whole metabolome with a single workflow.

Novel Aspect

Using a novel flow injection/SIM/ddMS² method on a hybrid quadrupole-Orbitrap instrument we developed a LC-MS² library, that can be used for Informed unTargeted Metabolomics.

TP-259 / LC-MS ANALYSIS OF NONRIBOSOMAL PEPTIDOME FROM A BRAZILIAN CYANOBACTERIUM NOSTOC SP. CENA543

Jouni Jokela⁽¹⁾ - *Lassi Heinilä*⁽¹⁾ - *Tânia Keiko Shishido*⁽¹⁾ - *Matti Wahlsten*⁽¹⁾ - *David P. Fewer*⁽¹⁾ - *Marli Fátima Fiore*⁽²⁾ - *Esa Haapaniemi*⁽³⁾ - *Perttu Permi*^(3,4) - *Kaarina Sivonen*⁽¹⁾

University of Helsinki, Department of Microbiology, Helsinki, Finland (1) - University of São Paulo, Center for Nuclear Energy in Agriculture, Piracicaba, São Paulo, Brazil (2) - University of Jyväskylä, Department of Chemistry, Jyväskylä, Finland (3) - University of Jyväskylä, Department of Environmental Science, Jyväskylä, Finland (4)

Keywords: Nostoc, cyanobacterium, nodularin, pseudospumigin, namalide

Introduction

Nostoc is a cosmopolitan genus of cyanobacteria found in soils and in symbiotic associations and a prolific producer of bioactive natural products. Nodularin (NOD) is a cyclic nonribosomal pentapeptide, hepatotoxic protein phosphatase inhibitor and tumor promoter. Spumigins are small linear nonribosomal peptide protease inhibitors. Anabaenopeptins are cyclic nonribosomal hexapeptides. Namalides are cyclic tetrapeptides structurally as anabaenopeptin.

Methods

LC-MS analyses were performed with an Agilent 1100 Series LC/MSD Ion Trap XCT Plus System (Agilent Technologies, Palo Alto, CA, USA). High accuracy mass was measured by UPLC-ESI-QTOF mass spectrometer Synapt G2 Si HDMS (Waters Corp.).

Results

In this study we analyzed non-ribosomal peptidomes with LC-MS methods from 18 cyanobacterial genera that had been collected from Brazilian saline-alkaline lakes of the Nhecolândia, Pantanal wetland area. Our objective was to discover novel bioactive compounds. We found a cyanobacterium strain *Nostoc* sp. CENA543 which produced 18 different peptide variants from which 13 were novel. Peptides belong to four different peptide families namely nodularins [1], pseudospumigins [1], anabaenopeptins [2] and namalides [2]. Pseudospumigins constitute a novel trypsin inhibiting peptide family. Pseudospumigin A inhibited trypsin (IC50 4.5 μ M after 1 h). Hepatotoxic nodularin was produced at exceptional high amounts, at same levels as many *Nodularia spumigena* strains isolated from blooms from different geographical locations produces nodularin. Together our results demonstrate that the nonribosomal peptidome of *Nostoc* sp. CENA543 analyzed at gene and metabolite level is strikingly similar to *N. spumigena* but contains also clear differences.

Conclusions

Natural products are an important source of antimicrobial, antifungal, anticancer, immunosuppressant and other bioactive molecules. We demonstrate that *Nostoc* sp. CENA543 produces one phosphatase inhibitor peptide group and three protease inhibitor peptide group pseudospumigins, anabaenopeptins and namalides with 13 novel peptide structures.

Novel Aspect

Complete genome contains biosynthetic gene cluster only for anabaenopeptin so namalides are a module skipping products from the anabaenopeptin biosynthetic pathway.

References

1. Jokela J., Heinilä LM. P., Shishido T. K., Wahlsten M., Fewer D. P., Fiore M. F., Wang H., Haapaniemi E., Permi P., Sivonen K. *Front of Microbiol*, 8:1963 (2017).
2. Shishido T. K., Jokela J., Fewer D. P., Wahlsten M., Fiore M. F., Sivonen K. *ACS Chem Biol*, 12(11), 2746-2755 (2017).

TP-260 / COMPARATIVE EVALUATION OF DIFFERENT GINSENG BERRY SPECIES BASED ON UPLC-QTOF/MS AND HR-MAS NMR METABOLITE PROFILES

Dae Young Lee ⁽¹⁾ - ***Jae Won Lee*** ⁽¹⁾ - ***Bo Ram Choi*** ⁽¹⁾ - ***Dahye Yoon*** ⁽²⁾ - ***Young Chang Kim*** ⁽¹⁾ - ***Suhkmann Kim*** ⁽²⁾

National Institute of Horticultural and Herbal Science, Department of Herbal Crop Research, Eumseong, Republic of Korea (1) - Center for Proteome Biophysics and Chemistry Institute for Functional Materials, Pusan National University, Busan, Republic of Korea (2)

Keywords: ginseng berry, ginsenoside, UPLC-QTOF/MS, HR-MAS NMR

Introduction

Many studies have previously reported the pharmacological activities of ginsenosides from *Panax ginseng* C. A. Meyer. The interests in therapeutic potential of ginsenosides also have increased. Ginsenosides are distributed in several parts of *P. ginseng*, including not only root but also leaf and berry. Different parts of *P. ginseng* may contain distinct composition of ginsenosides and pharmacological activities [1]. Recent studies have indicated the utility of ginseng berry as a good ingredient for various ginsenosides [2,3]. Thus, in the herbal market and food industry, it is critical to evaluate and control the quality of ginseng berry. The profiling of ginsenosides and primary metabolites will be applicable for this purpose.

Methods

The ginseng berry samples were obtained from four different *P. ginseng* species (*i.e.* Sunone (SW),

Gumpoong (GP), Yunpoong (YP), Gopoong (GoP)). From the extract mixture, ginsenosides were separated by using UPLC system with ODS column. Exact mass measurements based on QTOF/MS were also effective for the reliable identification of individual ginsenosides. Also, different species samples of *P. ginseng* berry were analyzed through global metabolite profiling by HR-MAS NMR spectroscopy.

Results

In the negative ion mode of MS, ginsenosides were detected as $[M+COOH]^-$ ions. As a result, 35 ginsenosides were profiled from the four ginseng berries. The raw data of many metabolites including ginsenosides from four ginseng berries were subjected into the principal component analysis (PCA). In the PCA score plot, the four ginseng berries were well separated. The PCA loading plot also represented several ginsenosides as the major variables to discriminate the four ginseng berries. 16 primary metabolites of *P. ginseng* berry were identified by 1H NMR spectra. The PLS-DA clearly distinguished *P. ginseng* berry extract according to the four ginseng berries.

Conclusions

These results demonstrated that contents of ginsenosides differ depending on the species of ginseng berry, and UPLC-QTOF/MS and NMR based profiling method can be a good tool to evaluate the metabolites contents of ginseng berry. This method will help the quality control and effective use of ginseng berry in the food industry.

Novel Aspect

This study indicated that the proposed method based on the UPLC-QTOF/MS and NMR can be an effective tool for characteristics of Korean ginseng varieties.

References

1. Ahn I. O., Lee S. S., Lee J. H., Lee M. J., Jo B. G., J. Ginseng Res., 32, 15-18 (2008).
2. Attele A. S. et al. Diabetes, 51, 1851-1858 (2002).
3. Xie J. T. et al. Phytomedicine, 9, 254-258 (2002).

TP-261 / UNTARGETED METABOLOMICS WITH PARALLEL ACCUMULATION – SERIAL FRAGMENTATION (PASEF) ON A NOVEL TRAPPED ION MOBILITY MASS SPECTROMETER

Florian Meier⁽¹⁾ - **Catherine Vasilopoulou**⁽¹⁾ - **Karolina Sulek**⁽²⁾ - **Andreas-David Brunner**⁽¹⁾ - **Ulrike Schweiger Hufnagel**⁽³⁾ - **Sven Meyer**⁽³⁾ - **Aiko Barsch**⁽³⁾ - **Matthias Mann**⁽¹⁾

Max Planck Institute of Biochemistry, Proteomics and Signal Transduction, Martinsried, Germany (1) - NNF Center for Protein Research, University of Copenhagen, Copenhagen, Denmark (2) - Bruker Daltonik GmbH, Bremen, Bremen, Germany (3)

Keywords: metabolomics, lipidomics, ion mobility spectrometry, TIMS, PASEF

Introduction

The metabolome is a highly dynamic network of small molecules that represent not only the end-product of cellular metabolism, but also regulate and catalyze key biological processes. Here, we integrate Trapped ion mobility spectrometry (TIMS) and “Parallel Accumulation – Serial Fragmentation” (PASEF) [1, 2] to multiply the speed and sensitivity of MS2 acquisition in untargeted liquid chromatography – mass spectrometry-based lipidomics and metabolomics.

Methods

High-flow and nano-flow liquid chromatography was coupled online to a hybrid TIMS quadrupole time-of-flight mass spectrometer (Bruker timsTOF Pro) via electrospray ionization. Complex lipid and metabolite

extracts from tissue and plasma sample were measured in PASEF mode. Raw LC-TIMS-MS data were processed with MetaboScape (Bruker) using on four-dimensional feature recognition and ion mobility-enhanced compound identification from spectral libraries.

Results

IMS can separate even isobaric ions within ~100 ms, compatible with LC and TOF time scales. In TIMS this happens by the opposing forces of an electric field and the drag of a gas flow. In PASEF, the quadrupole switches rapidly within a single TIMS scan to capture as many precursors as possible for fragmentation without diminishing sensitivity. We found that in crude metabolite extracts, we could fragment up to 17 precursors in 100 ms TIMS scans, which increased the total number of data-dependent MS² scans by more than an order of magnitude. This extreme sequencing speed can be used in part to increase spectral quality for low-abundance precursors by re-sequencing. In terms of its physical principle, TIMS closely resembles conventional drift-tube IMS. This allows straightforward determination of collisional cross sections (CCS). We measured CCS values for a set of metabolites comprising vitamins, hormones, and lipids, which were in good agreement with literature values and highly reproducible with median coefficients of variation below 1%.

Conclusions

TIMS literally adds a novel dimension to conventional metabolomics LC-MS workflows that is only beginning to be explored. The speed of PASEF enables very high sample throughput. We anticipate that the accurate and precise measurement of CCS values in conjunction with highly sensitive PASEF acquisition will allow a deeper metabolomics and lipidomics characterization of, for example, clinical samples.

Novel Aspect

Robust high-speed and high-sensitivity metabolomics with online PASEF

References

1. Ridgeway, Park, Int. J. Mass Spectrometry, 425, 22 (2018).
2. Meier, Mann, Journal of Proteome Research, 14, 5378 (2015).

TP-262 / AUTOMATED PARALLEL DERIVATIZATION STRATEGY WITH BROAD METABOLITE COVERAGE COUPLED TO SWATH/MS DATA ACQUISITION FOR QUALITATIVE AND QUANTITATIVE ANALYSIS

David Ruskic⁽¹⁾ - *Maria Fernanda Cifuentes Girard*⁽¹⁾ - *Renzo Picenoni*⁽²⁾ - *Mario Francesco Mirabelli*⁽²⁾ - *Günter Böhm*⁽²⁾ - *Gérard Hopfgartner*⁽¹⁾

Life Sciences Mass Spectrometry, Department of analytical and Inorganic Chemistry, University of Geneva, Geneva, Switzerland (1) - CTC Analytics AG, Zwingen, Switzerland (2)

Keywords: Parallel derivatization, QUAL/QUANT, dansyl-derivatives, SWATH/MS, automation

Introduction

Diversity in metabolite polarity and electrospray ionization efficiency challenge the development of a single method. Chemical derivatization can significantly improve chromatographic retention time and MS response. We propose a workflow including an automated fast parallel derivatization of amines, phenols, aldehydes, alcohols and ketones followed by QUAL/QUANT SWATH/MS analysis for broad metabolite coverage.

Methods

Dansyl-chloride (DanCl), dansyl-hydrazine (Dan-N₂H₃) and their ¹³C label analogs were used as derivatization agents. Samples (49 analyte mix and urine) were derivatized with ¹²C reagents while standards were

derivatized with ^{13}C labeled reagents using PAL RTC autosampler (CTC Analytics). The two fractions were mixed and injected onto a column-switching LC system. MS acquisition was performed on a TTOF6600 (Sciex) using SWATH acquisition.

Results

Polar metabolites are converted to more hydrophobic products (DanCl and Dan- N_2H_3 derivatives), enabling them to be separated on reverse phase liquid chromatography. The presence of a basic p-amino group on both reagents increase the electrospray response factor by a factor of 10 to 200. Primary alcohols, phenols, primary and secondary amines are derivatized with DanCl and Dan- N_2H_3 transforms ketones and aldehydes to hydrazones. The automated parallel derivatization enables reproducible derivatizations in a single workflow including light labeling of sample and heavy ^{13}C labeling of metabolite mix in 15 minutes during the LC analysis of the previous sample. Collision induced dissociation generates fragments specific for the analyte, and for the light and heavy derivatized analytes and for tags (XIC $m/z=171$ and $m/z=173$). LC-SWATH/MS which collects all precursors and all fragments allowed the screening and the relative and absolute quantification (n=45) of metabolites in urine.

Conclusions

Labelled reagents enable to generate adequate standards for amines, phenols, aldehydes, alcohols and ketones. Automation was found to be key for reproducibility and is performed on-line prior LC analysis. In a batch, sample preparation and sample analysis are overlaid resulting in a significant gain of time. Dansyl-chloride and dansyl-hydrazine were found to be ideal for light/heavy labeling, improved LC retention, improved MS response and MS/MS tag

Novel Aspect

Automated parallel heavy || light derivatization of sample || standard mix using dansyl-chloride and dansyl-hydrazide for improved metabolites coverage with LC-SWATH/MS.

TP-263 / TARGETED METABOLOMICS APPLIED TO HEART FAILURE PATIENTS – A PILOT STUDY

Ralf Nilsson⁽¹⁾ - **Lars Löfgren**⁽¹⁾ - **Filippos Michopoulos**⁽²⁾ - **Maria Ryaboshapkina**⁽¹⁾ - **Pia Davidsson**⁽¹⁾ - **Cecilia Linde**⁽³⁾ - **Bengt Persson**⁽⁴⁾

Astrazeneca, Translational science, Mölndal, Sweden (1) - Astrazeneca, Oncology IMED Biotech Unit, Cambridge, United Kingdom (2) - Karolinska University Hospital,, Department of Cardiology, Stockholm, Sweden (3) - Danderyd Hospital, Cardiology, Stockholm, Sweden (4)

Keywords: Targeted metabolomics, heartfailure, reduced ejection fraction (HFrEF), preserved ejection fraction (HFpEF)

Introduction

Targeted metabolomics was performed at Astrazeneca on a pilot study with plasma samples from heartfailure (HF) patients in the 4D heartfailure PREFERS project for improvement of heartfailure management and research in Stockholm, initiated by the Karolinska Institute and Stockholm County. The patients were divided in three groups: HF with reduced ejection fraction (HFrEF, n=13), HF with preserved ejection fraction (HFpEF, n=22) and a HF control (n=41) group. All 76 patients had coronary artery bypass surgery performed.

Methods

The effective analysis of polar ionic metabolites by LC-MS, represents an analytical challenge for metabolic profiling. Ion-pair-chromatography (IPC) coupled to negative-electrospray-ionization tandem-mass spectrometry (1) was used to profile 100+ endogenous metabolic intermediates in human plasma. A

complementary method based on hydrophilic interaction chromatography (HILIC) positive-electrospray-ionization tandem-mass spectrometry was developed covering another 70+ metabolites. By combining these two platforms metabolites representing amino acids, phosphorylated sugars, citric acid cycle intermediates, glycolytic intermediates, acylcarnitines, lysophospholipids, and nucleosides were covered.

Results

There was no clear metabolite phenotype separation between normal, HFpEF and HFrEF patients based on targeted metabolomics but two metabolites were differentially abundant between HFrEF patients and controls and survived adjustment for multiple testing. There was also significant pair-wise separation between groups for more than 15 metabolites.

Conclusions

In conclusion, targeted metabolomics was successfully applied to the heart failure pilot study. results were in-line with published data, confirming the capability of the platform to support future work on an expanded cohort aiming at identifying dysregulated metabolic pathways to advance the understanding of heart failure pathophysiology in search for new treatments and biomarkers.

Novel Aspect

The application of combined targeted metabolomic platforms on samples from heart failure patients with preserved or reduced ejection fraction.

References

1. Michopoulos, Fet al. Targeted profiling of polar intracellular metabolites using ion-pair-high performance liquid chromatography coupled to tandem mass spectrometry: applications to serum, urine and tissue extracts. J Chromatogr. A. 2014 Jul 4, 1349, 60-8.

TP-264 / PHOSPHOLIPID REMOVAL BY SOLID PHASE EXTRACTION FOR UNTARGETED METABOLOMICS: IMPROVED SENSITIVITY FOR POLAR METABOLITES IN BRAIN

Päivi Pöhö⁽¹⁾ - Anu Vaikkinen⁽¹⁾ - Elena Scurti⁽¹⁾ - Mikko Airavaara⁽²⁾ - Tapio Kotiaho⁽¹⁾ - Risto Kostiainen⁽¹⁾

University of Helsinki, Division of Pharmaceutical Chemistry and Technology, Helsinki, Finland (1) - University of Helsinki, Institute of Biotechnology, Helsinki, Finland (2)

Keywords: Untargeted metabolomics, HILIC, solid phase extraction, phospholipid

Introduction

Untargeted UHPLC-ESI-HRMS metabolomics of polar brain tissue metabolites is challenging because brains contain large amounts of phospholipids. As phospholipids ionize easily, they cause signal suppression of polar metabolites and detector overloading at large injection volumes. Here we compare four solid phase extraction cartridges for phospholipid removal from brain extracts.

Methods

Rat brain proteins were precipitated with methanol, and supernatants subjected to phospholipid removal by pass through SPE. Non-encapped, octadecyl functionalized silica (C18), zirconia coated silica (Hybrid SPE), and two types of hydrophilic-lipophilic balance polymer cartridges (HLB and PRiME HLB) were compared for the lipid removal. The analytes were separated by HILIC with BEH Amide column, and detected in negative ion mode with ESI-Q-TOF-HRMS.

Results

Amount of phospholipids was significantly reduced with all sorbent materials, and subsequently, concentration and injection of larger volumes of sample was enabled. However, the removal was less

efficient with PRiME HLB compared to the other cartridges. The C18 and HLB sorbents removed efficiently also other lipid types from the samples, such as fatty acids, sulfatides, and ceramides. This is expected based on their non-polar retention mechanisms. PRiME and Hybrid SPE showed variable removal efficiency for the lipid classes that do not contain the phosphate group. In addition to lipids, some analytes of interest were retained in the sorbents. For example, in the case of HybridSPE, lactic acid, glucose, the dipeptide neurotransmitter *N*-acetylaspartylglutamic acid, and its precursor *N*-acetylaspartate were not recovered from the samples, while the C18 sorbent retained GABA completely. Extraction RSDs were mainly below 20 % for C18 and HLB, whereas PRiME and HybridSPE showed larger variance.

Conclusions

All tested sorbent materials are able to reduce the amount of phospholipids in methanolic brain extracts with very simple pass through protocols and thus improve sensitivity towards polar metabolites. However, some metabolites are lost in the process. The HLB sorbent provided best combination of lipid removal efficiency, repeatability, and metabolite recovery for brain metabolite analyses.

Novel Aspect

Systematic evaluation of solid phase extraction sorbents for phospholipid removal for the untargeted metabolomics of polar brain metabolites.

TP-265 / COFFEE METABOLOMICS BY LC-HRMS AND MOLECULAR NETWORKING

Claudia M. Rezende⁽¹⁾ - ***Najla M. D. Sá***⁽²⁾ - ***Rafael G. Costa***⁽³⁾

Federal University of Rio De Janeiro, Chemistry Institute, Rio De Janeiro, Brazil (1) - Ufrj, Chemistry Institute, Rio De Janeiro, Brazil (2) - Federal University of Rio De Janeiro, Ladetec-Chemistry Institute, Rio De Janeiro, Brazil (3)

Keywords: Metabolomics, coffee, mass spectrometry, molecular networking

Introduction

Coffee is one of the most appreciated beverages in the world. It has a unique flavor and aroma, resulted from a complex chemical composition of the beans and influenced by several factors [1]. The beverage quality is directly related to the chemical composition of the beans and studies about quality x chemical composition have been carried out for decades. However, they usually focus on a small number of metabolites, not covering the complex chemical composition of the beans.

Methods

Raw (unroasted) coffee beans with different beverage qualities were extracted comprehensively and exhaustively using a two-step solvent extraction and analyzed using a non-target method by liquid chromatography-data dependent high-resolution mass spectrometry (LC-HRMS ddMS2) in the Thermo QExactive Plus. The obtained files were converted to the mzXMLopen format and submitted to target and non-target metabolomic analysis as well as to molecular networking using the GNPS online platform [2].

Results

Sixty-four substances from different chemical classes including xanthines, chlorogenic acids, diterpenes, serotonin amides, among others. were identified using the *m/z* value with high accuracy (ppm error < 5ppm), fragmentation profile and isotopic pattern as the criterion for identification. The GNPS online platform allowed the creation of molecular networks containing molecular families according to the similarities in their chemical structure, helping in the better understanding of the chemical composition of the beans. By PCA and PLS-DA multivariate data analyses sucrose and the lipid fraction showed high correlation with beans of better quality, while substances in coffee wax were correlated to low quality coffees.

Conclusions

With this work it was possible to develop a comprehensive methodology based on metabolomics to study and better understand the chemical composition of coffee beans and its relation to the quality of the beverage. So, it may help coffee growers, cooperatives and other professionals in the coffee chain to produce better quality products.

Novel Aspect

We have developed a comprehensive metabolomic study identifying a large number of coffee compounds from the non-volatile fraction together with the creation of molecular networks based on tandem mass spectrometry data.

References

1. Wintgens, J. N. Coffee: Growing, Processing, Sustainable Production; Wiley-VCH: Weinheim, Germany, 2004.
2. Wang, M. et al. Nature Biotechnology, 34, 828-837, 2016.

TP-266 / UNTARGETED METABOLOMICS OF NATIVE TASMANIAN PLANTS: COMPARISON WITH THEIR MEDICINALLY-KNOWN RELATIVES

Karolina Skraskova - Richard Wilson - Jason Smith

University of Tasmania, School of Chemistry, Hobart, Australia

Keywords: High resolution mass spectrometry, Liquid chromatography, Multivariate statistics, Plant analysis, Chemical map

Introduction

Scientific knowledge of chemistry and medicinal properties of native Australian flora is sparse and only occasionally applied in practice. This is especially valid for Tasmania where the information on the Aboriginal use of herbal medicine was poorly documented. This project compares the chemistry of selected native Tasmanian plants (NTP) and their taxonomical medicinal “relatives” as a first step towards a chemical map of Tasmanian flora.

Methods

Seven plants from three genera were included in this study. Each genus was represented by one (or two) NTP species and one non-native species. Extracts of plants were prepared using ultrasound aided extraction and subsequently analyzed by UHPLC coupled to an LTQ-Orbitrap. High resolution LCMS data were mined using multivariate statistics and targeted metabolomics approaches to compare chemical differences / similarities between the plants from the same genera.

Results

Principal component analysis and agglomerative hierarchical clustering helped to visualize the chemical relations between the different plants. Interestingly, all native/non-native plant pairs (from the same genera) showed a high degree of chemical similarity. The greatest similarity was detected between the two *Epilobium* spp. By searching for compounds previously detected in their non-native counterparts, targeted data mining provided insight into the chemistry of the native plants. We determined, for example, that both *Epilobium* spp. were rich in polyphenolic compounds. Some of these compounds were previously found to have neuroprotective, anti-inflammatory, antioxidant, anti-proliferative, and other pharmacological properties [1]. This finding suggests that native Tasmanian *Epilobium* species may hold promise for medicinal uses.

Conclusions

We demonstrate that high resolution LC-MS combined with multivariate statistics is a powerful analytical approach to explore the chemistry of previously unmapped plant species. This is the first comparative study of native Tasmanian & non-native medicinal plants from the same genera.

Novel Aspect

Novel experimental design and analytical approach based on ultrasonic extraction, UHPLC & HRMS was used to create a proof-of-concept chemical map of native Tasmanian plants.

References

1. Schepetkin I. A., Ramstead A. G., Kirpotina L. N., Voyich J. M., Jutila M. A., Quinn, M. T. *Phytotherapy Research*, 30, 1287 (2016)

TP-267 / HYDROPHILIC INTERACTION LIQUID CHROMATOGRAPHY – FIELD ASYMMETRIC WAVEFORM ION MOBILITY SPECTROMETRY – MASS SPECTROMETRY AS A PLATFORM FOR UNTARGETED METABOLOMIC STUDIES

Katarzyna Szykula - Matthew Turner - Colin Creaser - James C. Reynolds

Centre for Analytical Science, Department of Chemistry, Loughborough University, Loughborough, United Kingdom

Keywords: Liquid chromatography mass spectrometry, field asymmetric waveform ion mobility spectrometry, untargeted metabolomics

Introduction

The challenge of metabolic profiling requires the use of analytical platforms that maximize the peak capacity. Previous studies have shown that the hyphenation of highly orthogonal methods such as hydrophilic interaction liquid chromatography (HILIC) with field asymmetric waveform ion mobility spectrometry (FAIMS) and mass spectrometry (MS) increased the number of peaks detected [1,2].

Methods

The developed workflow included one-step sample preparation by protein precipitation with acetonitrile which was followed by a 13 minutes chromatographic run on a HILIC column. The elute from the column was ionized by positive ion ESI and separated by FAIMS with the compensation field (CF) scanned from -0.5 to 3.5 Td (DF 240 Td) in 1 second at an MS scan rate of 20 scans per second, resulting in the acquisition of 19 MS across each CF scan, with one scan for re-initialization.

Results

The described workflow was applied to the analysis of fresh and aged (72 hours) urine samples. The method has shown good reproducibility for LC retention times (%RSD 0.13-0.87) and peak areas (%RSD 4.46-11.31). Principal component analysis (PCA) demonstrated a clear separation of fresh and aged urine at different CFs. Several biomarkers (up- and down-regulated) were found and identified for aging of urine by the presented method. The highest change in abundance was observed for m/z 121.072 (urea dimer) which increased in concentration with the age of the sample.

Conclusions

The combination of HILIC with fast scanning FAIMS and MS was applied to the analysis of urine samples. This rapid method including a simple sample preparation shows the potential for untargeted metabolomics and biomarker discovery.

Novel Aspect

For the first time, the hyphenation of FAIMS with LC-MS was applied to untargeted metabolomic studies of urine samples.

References

1. Arthur K. L., Turner M. A., Reynolds J. C., Creaser C. S., *Analytical Chemistry*, 89 (6), 3452-3459 (2017).
2. Wernish S., Afshinnia F., Rajendiran T., Pennathur S., *Analytical and Bioanalytical Chemistry*, 410 (12), 2865-2877 (2018).

TP-268 / AN EXAMPLE OF VALIDATION OF AN MS-BASED METHOD

Burcu Binici - Ilker Un - Taner Gokcen - Gokhan Bilsel - Ahmet C. Goren

TUBITAK National Metrology Institute, Organic Chemistry Laboratory, 41470, Kocaeli, Turkey

Keywords: Method Validation, Sulfadiazine, Enrofloxacin, Bovine Tissue

Introduction

The measurement has an important role in our life. Many scientists put effort on this issue and try to accomplish and demonstrate the reliability of the measurement. For this purpose, the method validation process was created to evaluate the performance of an analytical method in terms of reliability, selectivity and sensitivity. In this study, the validation of the method that was developed to determine the antibiotics in food matrix was examined as a case study.

Methods

Many parameters like linearity-working range, limit of detection (LOD), limit of quantification (LOQ), precision, trueness and selectivity were evaluated during the method validation procedure. As a part of quality of measurement, the uncertainty value was calculated based on the model equation approach and traceability was established by purity assessment of the calibrants. In analytical process traditional solvent extraction and LC/MS-MS was used.

Results

In method validation process, method gave linear response between 0.04-6.5 mg/kg mass fraction levels. LOD was obtained as 0.005 mg/kg and LOQ was obtained as 0.017 mg/kg for both compounds. Trueness was evaluated via recovery experiments and for high and low mass fraction level in working range. Recovery value was obtained as 98.11% for sulfadiazine and 99.96% for enrofloxacin. Precision was evaluated in-a-day and intr-a-day for high and low level mass fraction. RSD% was obtained between (2.07-3.86)% for sulfadiazine and (2.52-3.23)% for enrofloxacin. Traceability of the method was performed by purity assessment of calibrants with qNMR. Purity value was obtained as (99.993±0.194)% for sulfadiazine and (99.521±0.230)% for enrofloxacin. Uncertainty evaluation was performed by model equation (bottom-up) approach. Relative uncertainty value is 11.4% for sulfadiazine and 11.1% for enrofloxacin.

Conclusions

The fitness of purpose of the method was evaluated by applying method validation considering the many parameters and so the developed method obtained the proof of performance. The comparability of the method was established by uncertainty and traceability assessment. In the case study a method was developed to determine high-polarity analytes ($pKow \geq -2$), sulfadiazine and enrofloxacin, in high fat and high protein food matrix.

Novel Aspect

This study presents an example of full method validation applied on simple, selective and reliable analytical method to measure sulfadiazine and enrofloxacin in bovine tissue.

WPS-S01 - Native MS and structural proteomics

Chairs: Frank Sobott, Nathalie Norais

WP-1 / Acidic XPlex: a multiplex alternative to NHS esters for faster cross-linking reactions at lower temperatures

Bruno Cesar Amaral⁽¹⁾ - ***Fabio Cesar Gozzo***⁽¹⁾ - ***Diogo Borges Lima***⁽²⁾ - ***Paulo Costa Carvalho***⁽³⁾

Dalton Mass Spectrometry Laboratory, University of Campinas, Institute of Chemistry, Campinas (1) - Mass Spectrometry for Biology Unit, Institut Pasteur, Paris (2) - Group for Computational Mass Spectrometry & Proteomics, Carlos Chagas Institute, Fiocruz Paraná, Curitiba (3)

Keywords: Cross-linking, XPlex, Adipic Acid, Multiplex

Introduction

Traditionally, chemical cross-linking is performed using NHS esters that react with nucleophilic groups. However, they are easily hydrolysable and have low solubility (sulfonated versions are very expensive). Reaction time, temperature and coverage are also suboptimal. Here, we describe a new cross-linking chemistry (Acid XPlex) overcoming these limitations using cheap, easily available reagents.

Methods

Myoglobin and cytochrome C were subjected to cross-linking reaction using XPlex activated adipic acid at 5, 15 and 27 Celsius degrees, followed by tryptic digestion. Samples were subjected to MS and MSMS analysis a Q-TOF and an Orbitrap Fusion Lumos. Data were interpreted using Waters MassLynx and SIM-XL software. Cross-linking reactions using DSS were also performed for comparative reasons.

Results

After condition optimization steps, reaction of myoglobin using adipic acid at 27 Celsius degrees showed 9 cross-linker molecules incorporation per protein compared to 6 of DSS. This difference is even more pronounced at lower temperatures where at 5 degrees only 2 DSS incorporations were observed, compared to 8 of adipic acid. The same results were obtained when using cytochrome C as the target protein. The method also allowed for a reduced reaction time (from 2 hours using DSS to 30 minutes using adipic acid), Another advantage of XPlex chemistry is the broader reactivity as cross-links are observed between Lys/Lys, Lys/Ser, Ser/Ser, Lys/Glu, Lys/Asp, Ser/Glu and Ser/Asp residues, improving the protein surface coverage and number of distance constraints. Special algorithms have been implemented in SIM-XL software to efficiently identify these multiplex cross-links.

Conclusions

The Acid XPlex method presents several advantages over the use of NHS based cross-linkers, such as higher solubility, reduced reaction time, lower reaction temperatures, higher incorporation rates and, more importantly, broader reactivity and protein surface coverage, showing a potential to expand the use of XL-MS to proteins unstable in solution at higher temperatures or for prolonged periods.

Novel Aspect

Novel chemical cross-linking chemistry with improved incorporation rates and reactivity at low temperatures.

WP-2 / INTACT MASS ANALYSIS OF MACROMOLECULAR PEG-FAB CONJUGATE USING MULTI-MODE NATIVE LC COUPLED TO ULTRA-HIGH MASS RANGE NATIVE MS

Aaron Bailey⁽¹⁾ - **Guanghui Han**⁽²⁾ - **Jonathan Josephs**⁽¹⁾ - **Wendy Sandoval**⁽²⁾ - **Eugen Damoc**⁽³⁾

Thermo Fisher Scientific, Chromatography and Mass Spectrometry, San Jose (1) - Genentech, Inc., Protein Chemistry, South San Francisco (2) - Thermo Fisher Scientific, Chromatography and Mass Spectrometry, Bremen (3)

Keywords: Pegylation, PEG dendrimer, ion exchange, size exclusion, antibody

Introduction

Monoclonal antibody Fab domain fragments can be utilized therapeutically to bind to and sequester specific antigens to decrease activity of disease targets in the cells of patients. Large multi-pronged polyethylene glycol (PEG)-based structures can be utilized to 'bundle' together multiple Fab domains to yield a drug with highly concentrated therapeutic activity and increased bioavailability (1). PEG-Fab conjugates pose enormous potential for sample heterogeneity, and thus present a significant analytical challenge in upstream characterization. Native MS and other charge reduction strategies yield increased spectral separation to allow the accurate intact mass analysis of such highly complex spectra. We present a semi-automated workflow for analyzing high mass PEG-protein conjugates utilizing native ion exchange coupled directly to an ultra-high mass range Orbitrap MS system.

Methods

PEG-Fab conjugate was analyzed intact in native MS conditions using nanospray infusion, size exclusion chromatography (SEC), or ion exchange (IEC) chromatography. PEG-Fab conjugate was analyzed by SEC (Thermo Scientific MabPac SEC-1) or IEC (Thermo Scientific ProPac WCX) using isocratic conditions and pH gradient conditions, respectively. For SEC we used a mobile phase of 50 mM ammonium acetate pH 6.6. For IEC we used a gradient of pH 6.6 to 10.2 in a background of 50 mM ammonium acetate at a flow rate of 200 $\mu\text{L}/\text{min}$. LC-MS was performed using a Vanquish H-Class UHPLC with single wavelength UV detection coupled directly to an Orbitrap MS instrument capable of ultra-high mass range detection up to 80,000 m/z .

Results

A major aspect of characterizing PEG-Fab conjugates is determining the weighted average of Fab domains which are attached to an individual PEG core structure. In order to assess the degree of Fab-conjugation, PEG-Fab conjugate was analyzed by static nano-spray infusion using an ultra-high mass range Orbitrap MS system operated with nitrogen or xenon as the HCD trapping gas for cooling ions during transmission. PEG-Fab conjugate in native conditions resulted in a complex spectrum which did not yield sufficient definition in order to confidently determine accurate masses to measure the average ratio of Fab-to-PEG. Addition of the charge reducing reagent TMGN produced a more simplified spectrum spanning a range of 20,000 to 50,000 m/z which enabled deconvolution for accurate mass determination. These results showed an average load of approximately 7 Fab domains, which is consistent with a low resolution result attained by MALDI-TOF MS data. We coupled IEC directly to our Orbitrap MS system to allow charge-based separation and on-line native intact mass analysis. We found PEG-Fab conjugate isoforms elute as individual peaks in increasing order of Fab load. IEC-MS spectra became further simplified compared to direct nano-spray infusion. Isoform identity and abundance could be further confirmed and relatively compared by leveraging both chromatographic and mass spectral data for intact mass analysis.

Conclusions

Native SEC and native IEC coupled directly to UV and Orbitrap MS with ultra-high mass range detection provides a flexible platform for analyzing large (500 kDa) compounds which exists in very complex mixtures such as PEGylated species. SEC and IEC have unique strengths in chromatographically resolving complexity

which allow the Orbitrap MS system to obtain increasingly detailed information of molecular structure.

Novel Aspect

Intact mass analysis of a complex PEG-protein conjugate using native ion exchange chromatography ultra-high mass range Orbitrap MS

References

- 1 Zarzar J, Shatz W, Peer N, Taing R, McGarry B, Liu Y, Greene DG, Zarraga IE. *BiophysChem.* 2018 May;236:22-30

WP-3 / NATIVE MASS SPECTROMETRY ANALYSIS OF RECOMBINANT ANTIBODIES DIRECTLY FROM THE GROWTH MEDIUM

Gili Ben-Nissan⁽¹⁾ - ***Shay Vimer***⁽¹⁾ - ***Aliza Katz***⁽²⁾ - ***Shira Warszawsky***⁽¹⁾ - ***Ron Diskin***⁽²⁾ - ***Sarel Fleishman***⁽¹⁾ - ***Michal Sharon***⁽¹⁾

Weizmann Institute of Science, Department of Biomolecular Sciences, Rehovot (1) -- Weizmann Institute of Science, Department of Structural Biology, Rehovot (2)

Keywords: Native Mass Spectrometry, Crude Proteins Samples, Recombinant Antibodies

Introduction

Typically, native mass spectrometry (MS) requires costly and labor-intensive protein purification. Formerly we presented a novel approach for native MS analysis of overexpressed proteins directly from crude bacterial lysates [1]. Here we expand the method to the analysis of secreted recombinant antibodies directly from the growth medium of secreting HEK293 cells.

Methods

Growth media from adherent HEK293 cells grown in DMEM, or suspension-grown HEK293 cells grown in Freestyle medium were transfected with plasmids encoding light and heavy chains of various antibodies. Growth media containing the secreted antibodies were collected and buffer exchanged into MS compatible buffers. Low-expressing samples were concentrated. The antibodies were analyzed using either an Orbitrap [2] or a Synapt G2 HDMS QTOF-based platform [3].

Results

Initially, we tested an antibody secreted from adherent HEK293 cells. The growth medium was 2x buffer exchanged into 1M ammonium acetate and concentrated 10x before analysis. The antibody produced a well-resolved charge state series which corresponded in mass to the glycosylated form. PNGaseF treatment directly in the medium efficiently removed the *N*-linked glycosylations. We next examined growth media from transfected suspension-cultured cells and showed that the secreted antibodies could be easily measured directly from the growth medium without any concentration. Antibodies could also be reduced and deglycosylated in the growth medium, for detailed structural analysis. For deeper characterization we compared wild-type and engineered antibodies, designed for increased affinity and stability, and demonstrated that the engineered antibodies were expressed to much higher levels and showed higher stability. CCS analysis demonstrated that the wild-type and designed antibodies exhibited identical collision cross-sections.

Conclusions

We showed that direct native MS can be efficiently used for the structural analysis of secreted antibodies from crude growth media. Analysis of the antibodies directly from the growth medium allowed us to rapidly determine the identity, expression level, modification level and overall structure of the antibodies. Moreover, we showed that our method is not dependent on a specific MS platform but rather it is a general

approach that can be applied on various instruments.

Novel Aspect

Direct native MS can be applied to the structural analysis of secreted antibodies from crude growth media, and can overcome the need for purification, saving time and effort.

References

1. Gan J., Ben-Nissan G., Arkind G., Tarnavsky M., Trudeau D., Noda Garcia Lianet., Tawfik D., Sharon M. *Analytical Chemistry*, 89: 4398-4404 (2017)
2. Ben-Nissan G., Belov ME., Morgenstern D., Levin Y., Dym O., Arkind G., Lipson C., Makarov AA., Sharon M. *Analytical Chemistry*, 89:4708-4715 (2017)
3. Hernández H., Robinson, C. V. *Nature Protocols*, 2:715–726 (2007)

WP-4 / COUPLING ONLINE SEPARATIONS TO NATIVE MS FOR THE CHARACTERIZATION OF DESIGNED PROTEIN COMPLEXES

Florian Busch⁽¹⁾ - **Mengxuan Jia**⁽¹⁾ - **Zachary VanAernum**⁽¹⁾ - **Aniruddha Sahasrabuddhe**⁽¹⁾ - **Zibo Chen**⁽²⁾ - **Scott Boyken**⁽²⁾ - **David Baker**⁽²⁾ - **Vicki Wysocki**⁽¹⁾

The Ohio State University, Department of Chemistry and Biochemistry, Columbus (1) - University of Washington, Institute for Protein Design, Seattle (2)

Keywords: native MS, protein design, protein complexes, structure screening, native separation

Introduction

De novo designed protein complexes have significant potential to serve as future drugs and nanomaterials. Whereas it is possible to design, express, and purify a great number of proteins and protein complexes in a short time, methods for their subsequent structural characterization on a large scale are missing. Here, we demonstrate how online separation coupled to native mass spectrometry can be used to rapidly assess key properties of designed protein complexes.

Methods

Protein samples were designed, expressed, and purified in the Baker lab. Gel filtration and ion exchange chromatography were performed with an Ultimate3000 HPLC system (Dionex) coupled to a modified QExactive EMR Orbitrap instrument (Thermo). UniDec [1] and Intact Mass (Protein Metrics) [2] were used for data deconvolution and analysis.

Results

We determined the oligomeric states and ligand-binding properties of several hundred designed protein complexes in TBS buffer by online gel filtration coupled to high-resolution mass spectrometry. Heterodimeric coiled-coil structures, with characteristics matching those intended by design were further evaluated for interaction partner specificity based on their tendency to exchange/ form off-target binding during a denaturation-refolding procedure. To deal with the complexity of such mixtures, we used ion exchange chromatography coupled to mass spectrometry for analysis and processed the data with Intact Mass. All-ion fragmentation was applied to confirm the protein complex identification. With several small mixing experiments including 6-16 heterodimers in each set, we filtered out 20 orthogonal heterodimers. For this set of 20 hetero-dimers, we found negligible subunit exchange during a denaturation-refolding procedure, consistent with the high specificity of the underlying designed protein-protein interactions.

Conclusions

Gel filtration coupled to mass spectrometry can be used to rapidly perform the initial identification of protein designs suitable for further structural analysis. Complicated mixtures can be readily separated by ion exchange chromatography coupled to MS and identified by fragmentation. In all cases, deconvolution

software is highly beneficial to determine the identity and oligomeric state of the designed protein complexes with high-throughput.

Novel Aspect

On-line separation methods for native MS were developed to help with the determination of structural properties of designed proteins at high speed with minimal sample preparation requirements.

References

1. Marty MT., Baldwin AJ., Marklund EG., Hochberg GKA., Benesch JLP., Robinson CV., *Analytical Chemistry*, 87 (8), 4370–4376 (2015)
2. Bern M., Caval T., Kil YJ., Tang W., Becker C., Carlson E., Kletter D., Ilker Sen K., Galy N., Hagemans D., Franc V., Heck AJR., *Journal of Proteome Research*, 17 (3), 1216–1226 (2018)

WP-5 / STRUCTURAL CHARACTERISATION OF WILD-TYPE AND Δ N6 β 2-MICROGLOBULIN COMPARING HDX AND FPOP CHEMICAL LABELLING

Owen Cornwell - James Ault - Sheena Radford - Alison Ashcroft

University of Leeds, Faculty of Biological Sciences, Leeds

Keywords: FPOP, HDX, β 2 – microglobulin, protein, aggregation

Introduction

Fast Photochemical Oxidation of Proteins (FPOP) is a covalent labelling technique that utilises UV radiation to produce hydroxyl radicals which covalently label solvent-accessible side-chains. Following proteolysis, LC-MS/MS is used to localise and quantify the oxidation particular residues. Here, we assess the utility of FPOP and HDX as parallel techniques to probe protein structure and dynamics on Δ N6, an aggregation-prone variant of β 2-microglobulin.

Methods

HDX-MS experiments were carried out using the LEAP Technologies automated system in-line with an M class ACQUITY LC and HDX manager followed by a Synapt G2Si MS/MS operating in HDMSE mode. FPOP was performed employing a 248 nm CompexPro Excimer UV laser to generate hydroxyl radicals used to oxidise solvent accessible side-chains. This was followed by protein proteolysis and LC-MS/MS using a Q ExactivePlus Orbitrap.

Results

>95% sequence coverage was obtained for both proteins using either HDX or FPOP. In the regions directly surrounding the *N*-terminal truncation, HDX indicates significantly higher deuterium uptake for Δ N6 compared with wild-type protein at the earliest labelling time-point (30 seconds). Key differences in labelling by FPOP were similarly clustered around the *N*-terminal truncation of the protein.

Examination of FPOP LC-MS/MS data showed certain modified residues eluting at multiple retention times, corresponding to the various positional isomers expected from hydroxyl radical attack on aromatic side-chains. Quantifying these peaks separately, changes in the labelling of positional isomers between the two proteins were observed (consistent with published NMR structures), and considerable differences in labelling trends between individual positional isomers and the combined residue level data were evaluated.

Conclusions

The comparison of HDX and FPOP shows that they serve as informative parallel techniques to probe protein structure and dynamics. Moreover, FPOP can be used to interrogate structural changes at sub-residue level resolution. As protein aggregation is often caused by subtle structural changes including the exposure of hydrophobic side-chains (both of which are very sensitive to oxidative labelling), FPOP has proved invaluable to the study of protein aggregation.

Novel Aspect

A comparison between two different probes of protein solvent accessibility and dynamics offers greater understanding and sub-residue resolution.

WP-6 / DOUBLE MUTANT CYCLES IN THE GAS PHASE: MEASURING INTER-PROTEIN PAIRWISE INTERACTION ENERGIES FROM CRUDE CELL LYSATES BY NATIVE MS

Jelena Cveticanin⁽¹⁾ - *Ravit Netzer*⁽¹⁾ - *Galina Arkind*⁽¹⁾ - *Sarel Fleishman*⁽¹⁾ - *Amnon Horovitz*⁽²⁾ - *Michal Sharon*⁽¹⁾

Weizmann Institute of Science, Department of Biomolecular Sciences, Rehovot (1) - Weizmann Institute of Science, Department of Structural Biology, Rehovot (2)

Keywords: native MS, crude cell lysates, coupling energies, inter-protein hydrogen bonds

Introduction

A powerful method to determine the energetic coupling between amino acids is double-mutant cycle analysis (DMC) [1]. Here, we show that double-mutant cycle analysis by native mass spectrometry [2] can be carried out for interactions in crude *Escherichia coli* cell extracts, thereby providing estimates for coupling energies under crowding conditions and obviating the need for protein purification and generating binding isotherms.

Methods

The applicability of this native MS-based approach was tested for three different double-mutant cycles using the 1:1 complex between E2 colicin endonuclease (colE) and the Im2 immunity protein (Im). All proteins comprising the cycle, the two wild-type proteins and the two mutant proteins were co-expressed in *E. coli* using the pRSFDuet expression plasmid. Native MS experiments were performed using an Orbitrap EMR platform.

Results

The charge series of the individual complexes forming each of the cycles could be well resolved in the crude lysates due to the overproduction of these proteins in comparison to the endogenous bacterial proteins and because of reduced binding of adducts. The relative concentrations in a cell lysate of the four co-existing complexes comprising each cycle were determined simultaneously from a single mass spectrum. Determining the concentrations of the complexes (without the need to define the concentrations of the individual proteins) allowed us to calculate the coupling energies under crowding conditions. We also determined the coupling energies using purified proteins in a buffer solution as a control. The values of the coupling energies determined for all three pairwise interactions in lysates are found to be very similar to the one determined for these interactions using purified proteins in buffer, thereby indicating that crowding has little effect on the strength of inter-molecular hydrogen bonds.

Conclusions

Double-mutant cycle analysis by native mass spectrometry can be carried out for interactions in crude *Escherichia coli* cell extracts, thereby providing estimates for coupling energies under crowding conditions and obviating the need for protein purification and generating binding isotherms. Our results indicate that inter-molecular hydrogen bond strengths are not affected by crowding conditions.

Novel Aspect

Our approach enables to characterize pairwise interactions by native MS of crude cell lysates obviating the need for any protein purification.

References

1. Horovitz A., Double-mutant cycles: a powerful tool for analyzing protein structure and function. *Fold &*

Des, 1 (1996).

2. Sokolovski M., Cveticanin J., Hayoun D., Korobko I., Sharon M., Horovitz A., Measuring inter-protein pairwise interaction energies from a single native mass spectrum by double-mutant cycle analysis. Nat. Commun. 8, 212 (2017).

WP-7 / ADVANCES IN ORBITRAP™ INSTRUMENTATION FOR NATIVE TOP-DOWN ANALYSIS OF NON-COVALENT PROTEIN COMPLEXES

Eugen Damoc⁽¹⁾ - ***Rosa Viner***⁽²⁾ - ***Albert Konijnenberg***⁽³⁾ - ***Kyle Fort***⁽¹⁾ - ***Maria Reinhardt-Szyba***⁽¹⁾ - ***Mikhail Belov***⁽¹⁾ - ***Alexander Makarov***⁽¹⁾

Thermo Fisher Scientific, Bremen, Bremen (1) - Thermo Fisher Scientific, San Jose, San Jose (2) - Thermo Fisher Scientific, Eindhoven, Eindhoven (3)

Keywords: Native MS, native top-down, pseudo-MS3, Orbitrap, Ultra-High Mass Range

Introduction

Native top-down studies of intact protein complexes have been reported since the early 1990's, but their characterization using MS3 have only recently been reported and most work has been done on homomeric assemblies. In this work we explore new ways for extending native top-down performance to allow interrogation of heteromeric protein assemblies like proteasome by top-down pseudo-MS3.

Methods

Experiments were performed using a modified Thermo Scientific™ Q Exactive™ Plus MS in which the transmission of high m/z ions and the ability to perform pseudo-MS3 scans for native top-down analysis were improved by implementing several hardware and software modifications (see the results section). Pyruvate Kinase, GroEL, rabbit 20S proteasome, LmrP membrane protein and AmtB membrane protein complex were used as model systems.

Results

Advances towards native MS and native top-down analysis of heteromeric protein complexes became possible only after addressing several major technical challenges. First, implementation of 'in-source trapping' capability addressed the insufficient or poorly controllable desolvation issue and allowed significantly improved fragmentation into subunits and stripped complexes in the inject flatapole region. Second, the reduction of RF frequencies on all RF guides and the mass filter was aimed at increasing the transmission of high m/z ions. Third, the adjustment of the voltage ramp rate on the central Orbitrap electrode facilitated successfully transmitting the high m/z ions from the C-trap into the Orbitrap analyzer. The improvements in performance afforded by these modifications were demonstrated in a series of pseudo-MS3 experiments performed for several homomeric and heteromeric intact soluble protein complexes and membrane protein complexes.

Novel Aspect

Improved native top-down analysis of non-covalent protein complexes

WP-8 / ARCHITECTURE OF EUKARYOTIC MRNA 3'-END PROCESSING MACHINERY - THE ROLE OF CROSS-LINKING AND MASS SPECTROMETRY

Gianluca Degliesposti⁽¹⁾ - ***Mark Skehel***⁽¹⁾ - ***Ana Casanal***⁽¹⁾ - ***Ananthanarayanan Kumar***⁽¹⁾ - ***Chris Hill***⁽¹⁾ - ***Ashley Easter***⁽¹⁾ - ***Paul Emsley***⁽¹⁾ - ***Yuliya Gordiyenko***⁽¹⁾ - ***Balaji Santhanam***⁽¹⁾ - ***Jana Wolf***⁽¹⁾ - ***Katrin Wiederhold***⁽¹⁾ - ***Gillian Dornan***⁽¹⁾ - ***Carol V. Robinson***⁽²⁾ - ***Lory A Passmore***⁽¹⁾

Laboratory of Molecular Biology (MRC LMB), Laboratory of Molecular Biology (MRC LMB), Cambridge (1) - Chemistry Research Laboratory, University of Oxford, Oxford (2)

Keywords: XL-MS, Cryo-EM, pre-mRNA processing, structural proteomics

Introduction

Eukaryotic pre-mRNAs are processed at their 3' ends by the Cleavage and Polyadenylation Factor (CPF). CPF cleaves pre-mRNAs, adds a poly(A)tail and triggers transcription termination. It is unclear how the complex is assembled and coordinated. A Cryo-EM investigation of CPF aided by nanoESI-MS and cross-linking-MS (XL-MS) shines new light on its assembly and activity providing an example of how XL-MS complements other structural approaches such as Cryo-EM. 1-4

Methods

XL-MS experiments were performed on purified CPF complex, which is composed of four or five proteins: Cft1, Yth1, Fip1, Pfs2, with and without Pap1. The isotopically-coded NHS ester BS3 was used as cross-linker. Trypsin digests were fractionated at the peptide-level by SEC and fractions analysed by nano-scale LC-MS/MS on an Orbitrap Velos. Cross-links were identified using xQuest and spectra were manually validated.

Results

NanoESI-MS of CPF complex revealed an interaction network organized into three modules: nuclease, phosphatase and polymerase. The Poly(A) polymerase module of CPF was expressed in insect cells. The ~200-kDa recombinant four-subunit complex (Cft1, Pfs2, Yth1, Fip1) was imaged using cryo-EM at a resolution of 3.5 Å. Cryo-EM maps display only three ordered subunits (Cft1, Pfs2 and Yth1). AXL-MS investigation was performed to validate the structural model and to determine where Pap1 and Fip1 bind. Cross-links map agree with the atomic models, and the crystal structure of Pap1. Fip1 cross-links to the C-terminal part of Yth1 and the polymerase domain of Pap1. Pap1 also cross-links the C-terminal helical domain of Cft1, ZnF1 of Yth1, and the C-terminal region of Pfs2. Together, these data suggest that the flexible C-terminal half of Yth1 binds the intrinsically disordered protein Fip1, which in turn flexibly tethers Pap1 to the complex, allowing conformational freedom to add long poly(A) tails onto diverse RNA substrates.

Conclusions

Nuclease, Poly(A) polymerase and phosphatase activities of yeast CPF are organized into three modules. A partial cryo-EM structure of polymerase module was solved. XL-MS provided insights on the interactions of the whole polymerase module. These data, combined with in-vitro reconstitution experiments, show that the polymerase module brings together factors required for specific and efficient polyadenylation, to help coordinate mRNA 3'-end processing.

Novel Aspect

Combination of nanoESI-MS, cryo-EM, XL-MS and in-vitro reconstruction experiments described the assembly and activity of the eukaryotic mRNA 3'-end processing machinery.

References

1. Rivera-Calzada, A. Degliesposti, G et al. The Structure of the R2TP Complex Defines a Platform for Recruiting Diverse Client Proteins to the HSP90 Molecular Chaperone System. *Structure* 25, 1145–1152 (2017)
2. Holzer, S. Degliesposti, G et al. Crystal structure of the N-terminal domain of human Timeless and its interaction with Tipin. *Nucleic Acids Res.* 45, 5555–5563 (2017)
3. Fiedorczuk, K. Degliesposti, G et al. Atomic structure of the entire mammalian mitochondrial complex I. *Nature* 538, 406–410 (2016)
4. Casañal, A. Degliesposti, G et al. Architecture of eukaryotic mRNA 3'-end processing machinery. *Science* 358, 1056–1059 (2017)

WP-9 / THE INFLUENCE OF SUBSTRATES AND NUCLEOTIDES ON COMPLEX FORMATION IN ABC TRANSPORTERS – INSIGHTS FROM MASS SPECTROMETRY

Francesco Fiorentino - Jani Reddy Bolla - Shahid Mehmood - Carol V. Robinson

University of Oxford, Department of Chemistry, Oxford

Keywords: ABC importers; native mass spectrometry; ModBC; BtuCD

Introduction

ABC importers are membrane proteins that couple the hydrolysis of ATP to the translocation of different molecules across cellular membranes. Depending on the fold of the transmembrane domains, they are classified as type I or type II importers[1]. To clarify the mechanism of these proteins, some questions need to be answered, such as the influence of substrate and ATP on the interactions between the transporters and their cognate substrate-binding protein.

Methods

We selected two very well characterised ABC importers: the molybdate transporter ModBC-A and B12 transporter BtuCD as representatives of type I and type II importer families, respectively. We employed native mass spectrometry-based methods (MS), enabling us the study of membrane proteins and their ligand-binding properties [2].

Results

Following purification of BtuCD and BtuF, we optimised native mass spectrometry conditions and obtained mass spectra of the two proteins in detergent solubilised form. When a solution of BtuF was added to BtuCD, the recorded mass spectrum showed the formation of the full BtuCD-F complex. After addition of vitamin B12, the spectrum still displayed the presence of the entire BtuCD-F complex but also revealed an increase in the intensity of the population of the uncomplexed BtuCD. Conversely, the addition of molybdate led to ModBC-A complex formation, while no complex was observed in the absence of substrate.

Mass spectra of BtuCD-F recorded after ATP addition showed the presence of the full complex bound to two ADP molecules. Interestingly, we detected an increase in the intensity of ADP-bound BtuCD. As for ModBC-A, following the addition of ATP we observed ADP-bound species, but we did not detect any additional charge states corresponding to ModBC.

Conclusions

Our MS data indicate that both vitamin B12 and ATP disrupt the interactions between BtuCD and BtuF. In contrast, molybdate is needed to let ModA dock onto ModBC and ATP does not affect ModBC-A complex formation.

These observations provide critical insights into the mechanism of transport of these importers and allow us to propose new transport models that complement the previous ones and will serve as a prototype for more complicated systems.

Novel Aspect

For the first time, native mass spectrometry has been applied to unravel the mechanism of action of transporters and provides further insights into their transport mechanism.

References

1. Locher, K. P. et al., Nat. Struct. Mol. Biol. 23, 487-493 (2016).
2. Mehmood, S. et al., Annu. Rev. Phys. Chem. 66, 453-74 (2015).

WP-10 / CHEMICO-STRUCTURAL AND FUNCTIONAL CHARACTERIZATION OF FILGRASTIM ORIGINATOR AND THREE OF ITS BIOSIMILAR DRUGS.

Alessandra Gianoncelli - Michela Bertuzzi - Michela Guarienti - Maurizio Memo

University of Brescia, Department of Molecular and Translational Medicine, Viale Europa, 11 - 25123 Brescia - Italy

Keywords: Originator, biosimilar, MALDI TOF/TOF mass spectrometry, liquid chromatography, zebrafish

Introduction

The approval authorization of a biosimilar is based on strong comparability studies with its biological references product due to the complexity of the structure and of the nature of manufacturing process of biological drugs. The proof of similarities needs both the development of analytical workflow, to characterize the protein chemically, and development of clinical trials to evaluate the efficacy and the safety of drug.

Methods

Qualitative and quantitative analyses were performed by tandem time of flight matrix-assisted laser desorption/ionization mass spectrometry (MALDI-TOF/TOF-MS) and reverse phase liquid chromatography coupled with an UV-detector (RP-HPLCUV). The evaluation of the biological activity of the biodrugs, in particular, their effect on immune system cell proliferation, was studied in zebrafish embryo, as animal model, using an innovative functional assay.

Results

The comparison studies between the recombinant filgrastim and its biosimilars were performed RP-HPLC and MALDI-TOF-MS both on intact proteins and on digested proteins. In particular, analyses were carried out to evaluate the integrity of proteins, the absence of drug-related impurities, the similarity of amino acid sequence and the protein concentration. While qualitative analyses allowed us to establish the molecular weight of the biodrugs and to verify the identity of primary amino acid sequence, quantitative analyses confirmed the protein concentration in the pharmaceutical formulations. After chemical characterization of the biodrugs, the biosimilarity was assessed by in-vivo study, using zebrafish embryos as model, to verify the biological activity. Functional studies showed that all drug had similar capacity to activate the innate immune response in zebrafish. Overall, the obtained results allowed us to demonstrate both the structural and biological similarity of the biosimilar drugs in comparison with their originator.

Conclusions

The analytical analyses performed in this study highlighted the structural similarity of biosimilar drugs Nivestim™, Tevagrastim® and Zarzio® and the originator Granulokine®. This result was reinforced by a similar in vivo biological activity.

Novel Aspect

This work is focused on the comparison of the originator of filgrastim with its three biosimilars, to evaluate the structural similarity and biological activity.

References

1. Gianoncelli et al, 2018. Submitted

WP-11 / CHARACTERIZATION OF CALCIUM-MEDIATED ALPHA-SYNUCLEIN AGGREGATION PATHWAY AT THE MOLECULAR LEVEL

Jong Yoon Han - Hugh I. Kim

Korea University, Department of Chemistry, Seoul, Republic of Korea

Keywords: α -synuclein, Intrinsically disordered protein structures, amyloid fibrillation, ion-mobility mass spectrometry (IM-MS), solution small-angle X-ray scattering (SAXS)

Introduction

Human α -synuclein (α Syn) is an intrinsically disordered protein, whose aggregation is associated with the pathogenesis of several neurodegenerative diseases (e.g., Parkinson's disease). Although α Syn aggregation process in the patients' brain is not fully understood, dyshomeostasis of Ca^{2+} was suggested to induce the aggregation. In this study, we examined the detailed structural change of α Syn and proposed the mechanism of Ca^{2+} -mediated α Syn aggregation [1].

Methods

Ion mobility-mass spectrometry (IM-MS) and solution small-angle X-ray scattering (SAXS) were utilized for structural characterization of monomeric α Syn in the presence of Ca^{2+} . Molecular dynamics (MD) simulations were performed to generate structures that match the experimental data. Multiple biophysical techniques were also utilized for revealing aggregation kinetics and morphological properties of aggregates.

Results

First, we observed that Ca^{2+} accelerated the α Syn fibrillation. To understand the role of Ca^{2+} , we investigated the structure of monomeric α Syn in the presence of Ca^{2+} using solution SAXS and IM-MS. Our IM-MS and SAXS results imply that population of extended conformers increases when Ca^{2+} ions are bound to α Syn. Based on the representative structures generated using MD simulation and ensemble optimization method (EOM), we suggested that the 1) Ca^{2+} induces α Syn to form extended conformer by electrostatic interaction, and 2) the extended conformers are prone to aggregation because hydrophobic NAC domain is more exposed to water than compact conformers. In addition, we observed that fibrils attached to each other in the presence of Ca^{2+} . Our detailed examination showed that Ca^{2+} binds to acidic-rich C-terminal domains originated from different fibrils, provoking interfibrillar aggregation.

Conclusions

In this study, we characterized the Ca^{2+} -mediated aggregation pathway of α Syn. Our results showed that binding of multiple Ca^{2+} ions induces the structural change of monomeric α Syn to expose NAC domain which is crucial for fibrillation, and provokes interfibrillar aggregation via electrostatic interaction.

Novel Aspect

Because of the biological importance of Ca^{2+} in α Syn aggregation, our study would be valuable for understanding pathology of α Syn-related neurodegenerative diseases.

References

1. J. Y. Han, T. S. Choi., Hugh I. Kim., Sci. Rep., 8, 1895 (2018).

WP-12 / CHARACTERIZATION OF INTACT MONOCLONAL ANTIBODIES UNDER NATIVE AND REVERSE PHASE CONDITIONS USING HIGH RESOLUTION MASS SPECTROMETRY

Chris Nortcliffe ⁽¹⁾ - Ryo Yokoyama ⁽²⁾ - Sibylle Heidelberger ⁽¹⁾ - Kerstin Pohl ⁽³⁾ - Jörg Dojahn ⁽³⁾ - Ferran Sanchez ⁽⁴⁾ - Annu Uppal ⁽⁵⁾

SCIEX, Biologics, Warrington (1) - SCIEX, Biologics, Tokyo (2) - SCIEX, Biologics, Darmstadt (3) - SCIEX, Biologics, Barcelona (4) - SCIEX, Biologics, Haryana (5)

Keywords: Native, NIST, Trastuzumab, TripleTOF 6600

Introduction

A major requirement in the development of biopharmaceuticals is extensive characterization of the molecule at the intact level to ensure safety and efficacy. More recently characterization of monoclonal antibodies under native conditions is gaining importance over reverse phase conditions as it helps preserve the noncovalent interactions and retain the folded structure.

Methods

Two monoclonal antibodies were studied using SCIEX TripleTOF® 6600 coupled to high-flow liquid chromatography (LC). Reverse phase chromatography was performed with water and acetonitrile containing 0.1% formic acid gradient. For native analysis, isocratic gradient of 50mM ammonium acetate buffer was used. MS detection was carried out using high resolution TOFMS mode.

Results

Two antibodies Trastuzumab, and NIST standard were analyzed under native and reverse phase conditions. The TripleTOF® 6600 system yielded high resolution accurate mass data with excellent sensitivity and reproducibility for both the molecules analyzed. The traditional reversed phase chromatography involves separation based on organic solvent gradient resulting in denaturing conditions and a broad distribution of highly charged species at lower m/z ranges. Relatively lower charge states are observed with narrow distributions and higher m/z under native conditions thus resulting in higher spatial resolution and improved baseline resolved peaks.

Conclusions

The high quality spectra obtained with TripleTOF® 6600 system along with the deconvolution using BioPharmaView™ software enabled the accurate molecular weight determination of all the possible glycoforms with excellent reproducibility.

Novel Aspect

Evaluation of the TripleTOF® 6600 systems for the Intact Mass analysis of mAb's under native conditions.

WP-13 / NATIVE MASS SPECTROMETRY OF BIOMOLECULES: A COMPARISON OF NESI-MS AND LILBID-MS

Nils Hellwig - Oliver Peetz - Nina Morgner

Institute of Physical and Theoretical Chemistry, Goethe-University, Frankfurt

Keywords: Native MS, LILBID-MS, nESI-MS, Membrane Proteins, Non-covalent Complexes

Introduction

Native mass spectrometry has emerged as an important tool in structural biology and is routinely used to analyze soluble proteins. However the investigation of membrane proteins and their complexes remains challenging; until the beginning of 2017, only about 30 membrane proteins have been investigated, mostly using nESI-MS. In this area LILBID-MS has recently shown great potential, which is further examined here in a comparison of the capabilities of these methods.

Methods

We compared the performance of the complimentary methods nESI and LILBID in the analysis of the soluble protein Avidin (tetramer) and the dimeric membrane protein EmrE. In addition, we performed a comparative screening of different buffers, concentrations and additives to explore the current limitations of both MS techniques. The nESI measurements were performed on a Waters Synapt G2-S, while the LILBID

spectra were recorded on an in-house build reTOF spectrometer.

Results

LILBID and nESI show different dissociation pathways. Avidin tetramer dissociates via dimers in LILBID, while in nESI, a highly charged monomer exits the complex. For EmrE, the native dimeric state can be detected almost exclusively using LILBID, whereas in nESI the dimer appears only as a minor species.

While with nESI, in case of most membrane proteins, it wasn't possible to avoid a certain amount of dissociation, we were able to demonstrate a gentle analyte release via IR-laser desorption in LILBID-MS, preserving almost the entire intact protein complex.

LILBID demonstrates a slightly higher tolerance towards non-volatile buffers TRIS-HCl and HEPES. High buffer concentrations require higher energy input for desolvation, leading to partial dissociation of the protein complexes. Both methods show a significant peak broadening with increasing salt concentration, the detected complex mass also increases. [1]

Conclusions

Soluble proteins and protein complexes can be analyzed convincingly using both methods.

If more desolvation of the analyte is needed, as for non-volatile buffers and additives or detergents used to stabilize membrane proteins, LILBID has benefits due to its different ion release mechanism.

When the native state of membrane protein complexes is of interest, LILBID can reveal this routinely, whereas nESI has its advantages in areas, where a high resolution is needed.

Novel Aspect

This is the first study comparing the capabilities of the two complimentary native MS methods nESI and LILBID in the analysis of soluble and membrane protein complexes.

References

1. Peetz, O., Hellwig, N., Henrich, E., Mezhyrova, J., Dötsch, V., Bernhard, F., Morgner, N. JASMS, submitted

WP-14 / EVALUATION OF AN ISOTOPE-LABELED MS/MS-CLEAVABLE CROSS-LINKER FOR PROTEIN STRUCTURE ANALYSIS

Patrizia Springorum^(1,2) - *Michael Götze*⁽²⁾ - *Christoph Hage*⁽¹⁾ - *Christian Ihling*⁽¹⁾ - *Mathias Schäfer*⁽³⁾ - *Andrea Sinz*⁽¹⁾

Department of Pharmaceutical Chemistry and Bioanalytics, Charles Tanford Protein Center, Martin Luther University Halle-Wittenberg, D-06120 Halle (Saale), Germany (1) - Institute of Biochemistry, Charles Tanford Protein Center, Martin Luther University Halle-Wittenberg, D-06120 Halle (Saale), Germany (2) - Institute for Organic Chemistry, Department of Chemistry, University of Cologne, D-50939 Cologne, Germany (3)

Keywords: Cross-linking mass spectrometry, MS/MS cleavable cross-linker, isotope-labeled cross-linker

Introduction

Previous studies have shown the benefits of the amine-reactive, CID-MS/MS-cleavable cross-linker DSBU (disuccinimidyl dibutyl urea) for structural proteomics studies based on cross-linking/MS [1, 2]. To further facilitate the automation of cross-linking/MS experiments we now introduce a D12-stable isotope labeled version of the DSBU-linker combining the advantages of MS/MS-cleavable linkers and isotope labeling.

Methods

Non-labeled DSBU and 1:1-mixtures of DSBU-D0 and -D12 were used to cross-link bovine serum albumine as well as the E. coli ribosome. Cross-linked proteins were proteolyzed and analyzed by nano-HPLC/nano-ESI-Orbitrap MS/MS (higher energy collision-induced dissociation, HCD) on Orbitrap Fusion and Q-Exactive

Plusmass spectrometers using mass difference-triggered MS/MS experiments. MS and MS/MS data were automatically analyzed and annotated with MeroX [3].

Results

The basis for a reliable, targeted selection of cross-linked peptides for MS/MS experiments is a mass difference between two isotopic species, which is large enough (D12) to create non-overlapping isotope patterns of a single species as well as overlapping elution profiles of the D0 and D12 species during LC/MS analysis. We could show that the D12-labeled DSBU-linker meets both requirements. The observed average retention time differences between D0 and D12 species were only ca. 10 seconds when employing 90-minute LC gradients. This enabled the efficient selection of cross-linked products for fragmentation by HCD, compensating the dilution effect caused by the application of two different cross-linker species. Employing mixtures of deuterated and non-deuterated cross-linkers combined with mass difference-based, targeted fragmentation seems to be advantageous compared to standard data dependent experiments using only one cross-linker species (156 versus 116 cross-linking sites for BSA).

Conclusions

The application of mixtures of deuterated and non-deuterated DSBU-linker is advantageous for structural proteomics studies based on cross-linking/MS. In combination with the dedicated software tool MeroX [3], a basis for automated, robust, and reliable automated cross-linking/MS studies is provided.

Novel Aspect

A D12-stable isotope labeled version of the CID-MS/MS cleavable DSBU-linker allows for an efficient mass-difference based, targeted analysis of cross-linked products.

References

1. Mueller, M. Q. ; Dreiocker, F. ; Ihling, C. H. ; Schaefer, M. ; Sinz, A. *Analytical Chemistry* 2010, 82, 6958-6968.
2. Arlt C., Götze M., Ihling CH., Hage C., Schäfer M., Sinz A., *Analytical Chemistry* 2016, 88 , 7930–7937
3. Goetze, M. ; Pettelkau, J. ; Fritzsche, R. ; Ihling, C. H. ; Schaefer, M. ; Sinz, A. *Journal of the American Society for Mass Spectrometry* 2015, 26, 83-97.

WP-15 / DEEP BIOSIMILARITY ASSESSMENT BY MONITORING MULTIPLE CRITICAL QUALITY ATTRIBUTES OF AN INTACT MONOCLONAL ANTIBODY DRUG USING NATIVE IEC AND SEC COUPLED TO NATIVE ORBITRAP MS

Aaron Bailey⁽¹⁾ - *Terry Zhang*⁽¹⁾ - *Stephane Houel*⁽¹⁾ - *Shanhua Lin*⁽²⁾ - *Guanghui Han*⁽³⁾ - *Bobby Gamez*⁽¹⁾ - *Katie Southwick*⁽¹⁾ - *Wendy Sandoval*⁽³⁾ - *Jonathan Josephs*⁽¹⁾

Thermo Fisher Scientific, 355 River Oaks Parkway, San Jose, CA 95134, USA (1) - Thermo Fisher Scientific 1228 Titan Way | Sunnyvale, CA 94085 (2) - Department of Microchemistry, Proteomics and Lipidomics, Genentech Inc., South San Francisco, CA 94080, USA (3)

Keywords: antibody, ion exchange, charge variant, native MS, size

Introduction

Manufacturing of innovator biologics can be successfully mimicked to produce generic “biosimilar” drug products. In order to satisfy safety and efficacy requirements, a biosimilar drug must be reasonably comparable to an innovator. Comparability is directly assessed by measuring a panel of critical quality attributes (CQAs). Recently, a multi-attribute method has been demonstrated to measure several CQAs simultaneously using LC-MS peptide mapping data. Additional CQAs, such as charge variants or size variants are conventionally monitored at the intact level using ion exchange (IEC) or size exclusion chromatography

(SEC), respectively, coupled to UV detection. Native LC-MS approaches, such as SEC-MS and IEC-MS, can be inclusive with UV detection to combine the traditional gold standard chromatographic assays with accurate mass measurement to allow isoform-specific monitoring of multiple CQAs of intact therapeutic proteins.

Methods

Intact trastuzumab was separated by SEC or IEC directly coupled to MS. Mobile phases consisted of aqueous 50 mM ammonium acetate. SEC was performed using isocratic pH 6.6 mobile phase. IEC was performed using a pH gradient from 6.8 to 10.1. LC-MS was accomplished using a Vanquish H-Class UHPLC system with a variable wavelength detector directly connected to a Thermo Scientific™ Q Exactive™ HF-X Orbitrap™ mass spectrometer. Native LC-MS raw data were analyzed using a time-resolved deconvolution approach utilizing Sliding Window and ReSpect algorithms in Thermo Scientific™ BioPharma Finder™ software. Thermo Scientific Chromeleon software was used to integrate XIC peak areas and relatively compare CQAs of innovator and biosimilar drug samples.

Results

Native IEC-MS analysis of intact trastuzumab resulted in a chromatographic profile consistent with previously published data. Time-resolved deconvolution analysis of the innovator drug resulted in identification of several specific isoforms comprised of differential *N*-glycosylation in combination with a variety of low level charge-imparting PTMs, such as deamidation, C-terminal lysine removal, or sialylation of *N*-glycans. We confirmed that the main peak consisted of several isoforms corresponding to variable *N*-glycosylations of the trastuzumab amino acid sequence. When comparing the innovator and biosimilar products, we found marked differences in glycoform profile measured at the main peak. Furthermore we found measureable differences in several charge variant isoforms when comparing innovator and biosimilar. These differences observed were consistent with a comparison of these same two samples via MAM HRAM peptide mapping analysis. Native SEC-MS of intact trastuzumab innovator and biosimilar samples showed marked differences in the levels of fragmentation as well as aggregation in response to a forced degradation treatment.

Conclusions

Biosimilarity of antibody therapeutic drug products can be accessed at the intact level to assess size variants and charge variants using either SEC or IEC coupled to UV and MS detection. This method was effective and reproducible for determining detailed molecular differences between an innovator drug, Herceptin, and a biosimilar.

Novel Aspect

Method for monitoring multiple CQAs of an intact biotherapeutic compound to allow direct comparison of innovator and biosimilar products in a convenient and deeply insightful. Both SEC and IEC coupled to UV and native MS can be used

References

1. Harris, R. J., Kabakoff, B., Macchi, F. D., Shen, F. J., Kwong, M., Andya, J. D., Shire, S. J., Bjork, N., Totpal, K., Chen, A. B. *J Chromatogr B Biomed Sci Appl.* 752, 233-245 (2001)

WP-16 / STRUCTURAL COMPARISON AND EPITOPE MAPPING OF INNOVATOR AND BIOSIMILAR THERAPEUTIC ANTIBODIES

Cristina Lento - Kerene Brown - Derek Wilson

York University, Chemistry, Toronto

Keywords: native MS, hydrogen/deuterium exchange, protein therapeutics, microfluidics

Introduction

Vascular endothelial growth factor (VEGF) is recognized as an essential regulator of normal and abnormal blood vessel growth [1]. The humanized mAb bevacizumab (Avastin) targets VEGF and blocks angiogenesis, resulting in suppressed tumour growth [2,3]. In this study, native MS, IM-MS, and global, as well as bottom-up HDX experiments were conducted to determine similarities in the higher-order structures of various lots of biosimilars to its innovator drug.

Methods

The main isoform of human VEGF (165) was purchased from R&D Systems. Avastin and biosimilars were exchanged into 250 mM $C_2H_7NO_2$ (pH 7.4) for all native MS, IM-MS and global HDX experiments. VEGF and antibodies were exchanged into 150 mM $C_2H_7NO_2$ (pH 6.7) and HDX was performed on the VEGF/Antibody complex and VEGF only to determine the epitope sites. Time-resolved global and bottom-up HDX was conducted using an in-house rapid mixing device and microfluidic platform.

Results

Native MS revealed that all six biosimilar lots had a more heterogeneous pattern and an overall greater level of glycosylation compared to Avastin. IMS-MS revealed identical drift times for Avastin and the biosimilars, consistent with gas-phase structures having the same collision cross-section. In the global HDX analysis, no differences in deuterium uptake were observed between Avastin and the biosimilars between 0.046 – 3.77 seconds. To localize subtle conformational changes and elucidate epitope sites, local HDX was conducted on VEGF only and the VEGF/antibody complex. At 8 and 4 seconds, the significant decreases in uptake identified for VEGF+Avastin were located at a. a. 27-43, 91-103 and 105-119 containing most of its known epitopes. Similar decreases were observed for VEGF+Biosimilar at a. a. 91-103 and 106-122. Interestingly, at lower time points, VEGF+Biosimilar displayed significant increases in deuterium uptake at a. a. 19-22, 91-103 and 154-162, while VEGF+Avastin continued to see decreases at the epitope, albeit to a lower extent.

Conclusions

Although there were differences observed in the glycosylation patterns, this did not seem to have an effect on the overall solution-phase dynamics of the antibodies. Although the binding site is the same for Avastin and the biosimilar, differences in allosteric perturbations could be observed at earlier time points which may be the result of how the antibodies engage VEGF. These increases could be attributed to differences in the thermodynamics of binding.

Novel Aspect

Using native structural mass spectrometry techniques in tandem with time-resolved HDX for biosimilar characterization and epitope mapping.

References

1. Ferrara N., Hillan K. J., Gerber H. P., Novotny W., Nat Rev Drug Discov, 3, 391-400 (2004).
2. Hicklin D. J., Ellis L. M., J Clin Oncol, 23, 1011-1027 (2005).
3. Presta L. G., Chen H., O'Connor S. J., Chisholm V., Meng Y. G., Krummen L., Winkler M., Ferrara N., Cancer Res, 57, 4593-4599 (1997).

WP-17 / LILBID-MS REVEALS INFLUENCES OF THT AND COOMASSIE ON AGGREGATION KINETICS OF AMYLOID-BETA (1-42) AND α -SYNUCLEIN

Janosch Martin - Rene Zangl - Tobias Lieblein - Nina Morgner

Goethe University, Institute of Physical and Theoretical Chemistry, Frankfurt

Keywords: LILBID-MS, Thioflavin T, Coomassie, amyloid-beta, α -Synuclein

Introduction

Many proteins can spontaneously self-assemble into β -sheet-rich fibrillar aggregates, so called amyloids. Two examples are amyloid-beta (1-42) (A β 42) and α -Synuclein (α -Syn), the proteins of interest regarding Parkinson's disease (PD) and Alzheimer's disease (AD), respectively. [1] To analyze aggregation of amyloids, a commonly used method is the Thioflavin T (ThT)-assay. Coomassie - used for staining proteins - is another important dye for protein analysis.

Methods

In contrast to a ThT-assay laser-induced-liquid-bead-ion-desorption mass spectrometry (LILBID-MS) enables to take a detailed look at small oligomers. Furthermore there is no need to add analysis required substances. Due to this LILBID can be used to study and compare not only the oligomerization kinetic by time-resolved measurements, but also to reveal the effect commonly used dyes like ThT and Coomassie have on this aggregation kinetic.

Results

LILBID was employed to analyze and compare the aggregation kinetics of α -Syn and A β 42. Using LILBID it is possible to detect a range of α -Syn and A β 42 oligomers. Time-resolved measurements demonstrate a faster aggregation of A β 42 compared to α -Syn.

The fluorescence dye ThT which tags β -sheet structures of a protein is still under discussion due to the potential to affect the aggregation of some proteins, *e.g.* Amyloid-beta (1-40). [2] Therefore the influence of ThT on the aggregation behavior of α -Syn and A β 42 was investigated. It can be shown that ThT neither promotes nor inhibits the oligomerization of these proteins, as long as if it is purified.

Recent studies also indicate an influence of the protein staining dye Coomassie on α -Syn aggregation. [3] Due to this, LILBID was also used to study how Coomassie affects the amyloidogenic oligomerization of α -Syn and A β 42. We could show that Coomassie promotes the oligomerization in case of α -Syn, whereas A β 42-aggregation is not influenced.

Conclusions

Time-resolved experiments were performed to analyze the influence of ThT and Coomassie on α -Syn and A β 42. It could be shown that ThT does not influence the aggregation behavior of these proteins. In contrast, Coomassie promotes formation of off-pathway oligomers in case of α -Syn, whereas A β 42-aggregation is not affected. These are important factors, relevant for interpretation of results obtained with the respective analysis compounds.

Novel Aspect

Requiring no additional substances, LILBID enables to analyze the potential influence of fluorescence dyes commonly used for aggregation analysis (see ThT-Assay).

References

1. Hamley I. W., *Angew. Chem.*, 119, 8274-8295 (2007).
2. D'Amico M., Di Carlo M. G., Groenning M., Militello V., Vetri V., Leone M., *J. Phys. Chem. Lett.*, 3, 1596-1601 (2012).
3. Ahsan N., Siddique I. A., Gupta S., Surolia A., *European Journal of Medicinal Chemistry*, 143, 1174-1184 (2018).

WP-18 / QUANTIFICATION OF SUBSTRATE AND INHIBITOR BINDING TO WILD TYPE AND MUTANT ISOCITRATE DEHYDROGENASE (IDH1) BY NON-DENATURING MASS SPECTROMETRY

Victor Mikhailov - Shuang Liu - Christopher Schofield - James McCullagh

Department of Chemistry, University of Oxford, Oxford

Keywords: native MS, drug binding affinities

Introduction

A single point R132H mutation of human isocitrate dehydrogenase 1 (IDH1) is responsible for inhibition of epigenetic regulators and is associated with a variety of cancers including over 70% of gliomas and 30% of myeloid leukemias. Inhibitors of the R132H IDH1 isoform are currently in preclinical and clinical trials. We present 'native' (non-denaturing) MS data on the binding of several of these inhibitors to IDH1 and we compare the results with data from other methods.

Methods

An IDH1 / inhibitor solution in ammonium acetate buffer, at different enzyme/inhibitor ratios, was analysed using a Nanomate nESI source (Advion) coupled to a Synapt qTOF MS (Waters). Inhibitors tested included Novartis 224 and GSK 864, and AGI-5198, AG 120 and enasidenib (Agiros). The binding stoichiometries for the wild type (wt) and R132H IDH1 and the dissociation constant KD were defined and native MS data were compared with the data from UV-absorbance, ITC and NMR.

Results

Mass spectra were consistent with the dimeric state of IDH1, single-site binding for AGI-5198 and AG-120, and two binding sites for GSK864 and Novartis 224. Enasidenib showed almost no binding to R132H IDH1 and single-site binding to wild type IDH1. KD measurements for the enzyme/inhibitor complex ranged from the low μM for AG-120 and GSK 864 to high μM values for AGI-5198/wt IDH1. For AG-120, GSK864 and Novartis 224, the difference between affinities to the wild type and R132H mutant enzyme was small. A large difference in KD was found for AGI-5198 binding to wt and R132H IDH1: 580 μM and 80 μM , respectively. KD values for AG-120 and GSK864 were in the same low μM range as the KD values from our NMR measurements. KD values for AGI-5198 and Novartis 224 were much higher than those from NMR. This may be indicative of a difference between the inhibitors in the binding sites. Native mass spectra of the complexes between the two forms of IDH1 and their respective substrates were also recorded. Studies on the effect of inhibitors on these complexes are currently underway.

Conclusions

We have demonstrated the utility of 'native' mass spectrometry in quantification of drug binding to a human enzyme involved in the development of brain and other cancers. We have compared binding stoichiometries and affinities of various inhibitors, and now are evaluating the inhibitor effect on the enzyme/substrate binding. Most of the data corresponds well with the results from other methods, the rest can be rationalised in terms of different binding mechanisms.

Novel Aspect

Stoichiometry of binding of various inhibitors to the IDH1 mutant and wild type enzyme that is involved in brain cancer, and quantification of binding affinities.

WP-19 / AN OPTIMIZED ENRICHMENT STRATEGY FOR IMPROVED MASS SPECTROMETRY ANALYSIS OF CHEMICALLY CROSS-LINKED PEPTIDES

Rosa Viner⁽¹⁾ - *Leigh Foster*⁽²⁾ - *Erum Raja*⁽²⁾ - *Chris Etienne*⁽²⁾ - *Ryan Bomgarden*⁽²⁾ - *Madalina Oppermann*⁽³⁾

Thermo Fisher Scientific, CMD, San Jose (1) - Thermo Fisher Scientific, CMD, Rockford (2) - Thermo Fisher Scientific, CMD, Stockholm (3)

Keywords: Cross-linking, MS-cleavable crosslinkers, enrichment, protein-protein interactions

Introduction

Chemical cross-linking in combination with mass spectrometry is a powerful method to determine protein-protein interactions. However, this method suffers from low identification rates without enrichment/fractionation. In this study, we assessed azido MS-cleavable crosslinkers. Specifically, we evaluated different probes, capture resins, elution protocols and fractionation for increased cross-linked peptide identification by mass spectrometry.

Methods

Amine-reactive, MS-cleavable crosslinkers including DSSO and Azido-A-DSBSO[1] were used to crosslink standard proteins and cell lysates. For azido-containing crosslinkers, samples were labeled with different biotin probes before digestion and enrichment. Samples were analyzed by LC/MS on the Thermo Scientific™ Orbitrap™ Lumos™ mass spectrometer. Thermo Scientific™ Proteome Discoverer™2.2 software and XlinkX node were used for data analysis.

Results

A major challenge of crosslinking mass spectrometry is the relatively low abundance of crosslinked peptides after protein digestion. This results in very few or no crosslinked peptide identifications by MS if the sample complexity is greater than a few hundred proteins. Therefore, enrichment of crosslinked peptides is necessary to identify them in complex samples. To specifically enrich modified peptides, affinity groups such as biotin have been incorporated into the crosslinker structure. However, biotin-enrichment of crosslinked peptides is also challenging due to highly abundant mono-linked peptides, contamination by endogenous biotinylated proteins and/or biotin in elution buffers, and poor recovery from avidin/streptavidin resins. In this work we optimized protocol for an enrichable crosslinker, Azido-A-DSBSO. The total time of sample prep was reduced 3 times. To further separate monolinks from fully crosslinked peptides, we also developed a workflow using biotin enrichment followed by SCX fractionation.

Conclusions

This new XL-MS workflow allows effective protein cross-linking in vivo to capture authentic protein interactions, selective enrichment of cross-linked proteins and peptides to improve their detection. Our optimized enrichment crosslinked peptide workflow with fractionation increased our identification of crosslinked peptides by 20% compared to enrichment alone for BSA samples and 30% for complex cell lysate samples.

Novel Aspect

Optimized sample preparation and enrichment of MS-cleavable crosslinked peptides for mass spec analysis.

References

1. Kaake R. M., et al, Mol. Cell. Proteomics, 13(12), 3533-3543 (2014).

WP-20 / TOWARDS DEFINING OF THE MOLECULAR MATURATION PATHWAY OF THE MITOCHONDRIAL SULFHYDRYL OXIDASE ERV1

Kamila Pacholarz⁽¹⁾ - ***Xiaofan Tang***⁽²⁾ - ***Hui Lu***⁽²⁾ - ***Perdita Barran***⁽¹⁾

Manchester Institute of Biotechnology, The University of Manchester, Manchester (1) - School of Biological Sciences, The University of Manchester, Manchester (2)

Keywords: Enzyme maturation, protein oligomerisation, cofactor binding, Erv1, activated IM-MS

Introduction

Erv1 (Essential for Respiration and Vegetative growth) is a FAD-dependent sulfhydryl oxidase which

catalyses disulfide bond formation and is required for the import and folding of small cysteine-containing proteins in the mitochondrial intermembrane space (IMS) [1,2]. Our aim is to use hybrid MS to define the maturation mechanism: oxidative folding, cofactor-binding and oligomerisation; of this essential enzyme to benefit the design of better biocatalysts.

Methods

Experiments were performed on full length Erv1 and its isolated compact C-terminal domain. Native MS and collisional activation prior to ion mobility, (aIM-MS) were applied to provide insights into Erv1 oligomerisation and FAD-binding (Waters Synapt G2S). FAD was titrated into apoErv1 to study the assembly pathway. Charge reducing agents were used to mimic the solution like oligomeric order of Erv1.

Results

Native MS experiments have shown full length Erv1 (purified in presence of FAD) to exist in tetrameric, dimeric and monomeric state. Both, tetrameric and dimeric Erv1 display over bimodal charge state distribution, characteristic for ordered-disordered systems. Still, Erv1 displays unusually broad charge state distribution even for a disordered system [3]. Addition of charge reducing agent to the buffer resulted in shift of the oligomerisation equilibrium to narrow dimeric and tetrameric states and mimicked solution measurements. Highly conserved C-terminal domain of Erv1 displayed a narrow charge state envelope characteristic of folded proteins. aIM-MS shows how intrinsically unfolded Erv1 species are more prone to release of FAD prior to the subunit dissociation. The energy barrier for subunit dissociation in compact multimers is lower than the cofactor release barrier. Titration experiments of FAD into apoErv1, present mainly as apo monomer and somewhat as apo dimer, show the assembly into a fully occupied tetrameric Erv1 upon binding of the cofactor.

Conclusions

Findings from MS provide new insights to define the maturation mechanism of Erv1. Erv1 exists in several multimeric states in the presence and absence of FAD, with a preference towards monomeric species in absence of the cofactor. FAD binding is essential for coordination of monomeric and dimeric Erv1 into tetramers.

Novel Aspect

Application of native MS and top-down IM-MS to study enzyme maturation process, along with use of charge reducing agent to mimic solution conditions.

References

1. Lee, J. E. et al., FEBS Lett., 477, 62, (2000).
2. Guo, P. -C., J. of Bio. Chem., 287, 34961, (2012).
3. Beveridge, R. et. al., Anal. Chem., 86, 10979, (2014).

WP-21 / ENRICHMENT OF ELECTROCHEMICALLY CLEAVED PEPTIDES FOR MIDDLE-DOWN PROTEOMICS APPLICATIONS

Hjalmar Permentier - Tao Zhang - Xiaobo Tian - Rainer Bischoff

University of Groningen, Analytical Biochemistry, Groningen

Keywords: electrochemistry, protein digestion, spirolactone, biotinylation

Introduction

Enzymatic digestion of proteins is standard procedure in bottom-up proteomics. Alternatively, digestion can be achieved by electrochemical (EC) oxidation of proteins, leading to specific cleavage at the C-terminal side of Tyr and Trp residues. EC oxidation generates a mixture of cleavage and non-cleavage products.

Cleaved peptides have a unique spirolactone at the C-terminus, which can be employed for their enrichment [1,2].

Methods

Model peptides and small proteins were oxidized at 1-2 V in an electrochemical flow-through cell with a carbon or boron-doped diamond electrode. The cleavage products were specifically derivatized with amine-containing reagents, including amine-biotin, at different pH values in the presence of metal ions. Derivatized peptides were affinity-enriched and analyzed by LC-MS/MS.

Results

Model tripeptides including LWL were oxidized to optimize the cleavage conditions and the reaction conditions with amine-biotin. At basic pH the peptide-spirolactone reactivity is highest, and a hydrolysis side reaction can be prevented by stabilization with Cu(II) ions. Other metal ions do not confer stability. The optimized conditions were then applied to small proteins, including lysozyme, leading to extensive cleavage at most Tyr and Trp residues. Enrichment of the cleaved, labelled peptides with streptavidin beads allowed the identification of the proteins in the digest. The amine-biotin tag readily fragments giving rise to a peak at m/z 375 which was used to verify the presence of enriched peptides. Further improvement of the enrichment using direct coupling of peptide-spirolactones to activated resins is expected to enable analysis of more complex EC-cleaved protein mixtures.

Conclusions

Cleavage of proteins with multiple Tyr and Trp residues and efficient derivatization with amine-biotin in the presence of Cu(II) was achieved. Biotinylated peptides were enriched with affinity capture, and LC-MS/MS analysis was used to identify the protein. Specific biotin fragments enable verification of biotinylation.

Novel Aspect

An electrochemical protein digestion method was optimized for protein identification using peptide derivatization and enrichment.

References

1. Roeser J., Alting N. F. A., Permentier H. P., Bruins A. P., Bischoff R. P. H., Rapid Commun. Mass Spectrom. 27, 546 (2013).
2. Zhang T., Niu X., Yuan T., Tessari M., de Vries M. P., Permentier H. P., Bischoff R., Anal. Chem. 88, 6465 (2016).

WP-22 / Heterogeneous Protein Complexes In Synaptic Vesicles Mediate Signal Transduction In Neurons

Sabine Wittig - Caroline Haupt - Marcelo Ganzella - Susann Kostmann - Reinhard Jahn - Carla Schmidt

HALOmem, Charles Tanford Protein Center, Institute for Biochemistry and Biotechnology, Martin Luther University, Halle

Keywords: Cross-linking, native MS, protein networks, synaptic vesicles, SNARE complex

Introduction

Synaptic vesicles (SVs) are small storage organelles for neurotransmitters. They are densely packed with proteins and pass through a trafficking cycle including fusion with the presynaptic membrane [1]. Membrane fusion depends on formation of the SNARE complex which connects the two membranes [2]. We combine chemical cross-linking and native mass spectrometry (MS) to unravel protein assemblies in SVs as well as the interactions formed during membrane fusion.

Methods

SVs were purified from rat brain. Individual proteins were expressed in *E. coli* and purified via a His-tag. Proteins were cross-linked using BS3 cross-linker. After tryptic digestion, cross-linked peptides were enriched by gel filtration, analysed by LC-MS/MS on a Q Exactive Plus mass spectrometer and identified using pLink software [3]. Assembled protein complexes were further studied by native MS using a modified Q-ToF Ultima mass spectrometer [4].

Results

SVs were purified from rat brain and major proteins were identified. Cross-linking then delivered first interaction networks in intact SVs showing that Syb plays a central role in complex formation presumably due to its high flexibility. Using botulinum toxin to cleaved the soluble segment of Syb and thus removed it from the protein networks allowing us to identify additional protein interactions. To study Syb in detail, we purified different variants, including full-length Syb, the soluble domain and a shorter segment. Using cross-linking and native MS we found that all variants formed oligomers in solution; the shortest variant formed higher order oligomers suggesting unspecific aggregation. By incubating the soluble segments of SNARE proteins or the full-length proteins incorporated into liposomes (*i.e.* membrane mimics) we studied the intermediates of SNARE complex formation. Cross-linking and native MS revealed protein interactions and stoichiometries in these assemblies.

Conclusions

Our results provide novel insights into the heterogeneous protein interaction networks in SVs which are key to our understanding of signal transduction in the neuronal synapse. In addition, by employing individually purified proteins and variants thereof, we provide insights into the interactions formed in solution or when the proteins are incorporated into membrane mimics. We are thus able to monitor SNARE complex formation and uncover intermediate states.

Novel Aspect

We provide first insights into protein interactions in SVs and pave the way to study membrane fusion by combining complementary MS techniques.

References

1. Südhof TC., Nature, 375, 645-53 (1995).
2. Jahn R., Fasshauer D., Nature, 490, 201-7 (2012).
3. Yang B., Wu YJ., Zhu M. et al., Nat Methods, 9, 904-9 (2012).
4. Sobott F., Hernandez H., McCammon MG., Tito MA., Robinson CV., Anal Chem, 74, 1402-7 (2002).

WP-23 / HYDROXYL RADICAL FOOTPRINTING STUDY FOR STRUCTURAL CHARACTERIZATION OF ALUMINUM HYDROXIDE-ADSORBED CALMODULIN BY MASS SPECTROMETRY

Alessandro Vadi

GSK Vaccines S. r. l., Analytical Research and Development, Siena

Keywords: Structural MS, Calmodulin, Aluminum hydroxide, Adjuvant adsorption, Hydroxyl radical footprinting.

Introduction

Possible conformational modifications of antigens during aluminum hydroxide adjuvant adsorption are discussed in literature for formulated vaccine for its potential impact on antigen immunogenicity [1-4]. The development of a mass spectrometry (MS) based technique called hydroxyl radicals footprinting (HRF) for the direct characterization of adjuvant-adsorbed protein is discussed in this poster showing a case study

on calmodulin (CaM).

CaM is a small and well characterized protein that is part of calcium signaling pathways [5].

CaM is an acidic protein with two globular domains connected by a central linker. Each domain contains two EF-hand sites that bind intracellular Ca^{2+} ions. Upon Ca^{2+} binding the EF-hand domains undergo to a conformational rearrangement, exposing hydrophobic protein domains that can bind a wide range of target proteins.

Hydroxyl radical footprinting technique (HRF) [6] coupled to MS was chosen to detect conformational changes of CaM, if any, when absorbed on aluminum hydroxide and then exposed to Ca^{2+} ions. Hydroxyl radicals are highly reactive species that chemically modify the surface accessible amino acid side chains with a chemical susceptibility to oxidation (+ 16 Da mass shift, or +32 or +14 depending on the residue), enabling monitoring of a wide number of target sites and consequently providing structural information.

Methods

Sample treatment protocol was developed ex novo to be compatible with aluminum hydroxide adjuvant. Fenton chemistry catalytic reaction [6] was used to generate hydroxyl radicals in solution and regenerate the Fe(II) species with sodium ascorbate. The three Fenton chemistry reagents (hydrogen peroxide, Fe(II)-EDTA, and sodium ascorbate) were added together with the sample. Hydroxyl radicals reaction with CaM was tested at 4 different time points (0, 2, 5, and 10 minutes) then reaction was quenched with thiourea. Subsequently protein were desorbed from adjuvant to allow the enzymatic digestion. Protein desorption was obtained by overnight incubation of the samples with phosphate buffer. Tryptic peptides are then analyzed on Q-Exactive Plus mass spectrometer (Thermo Fisher Scientific) by LC-MS/MS DDA top 5 experiment.

Data were processed to enable relative quantification of the oxidized methionine-residue. The 9 methionine residues distributed along the protein sequence of CaM were selected as probes for structural modification probing. Protein conformational characterization was correlated to the different solvent exposure of the methionine residues calculated on their oxidation amount linear regression slopes at different time points.

Results and Conclusions

Nine methionine residues (belonging to the protein sequence) were monitored on 4 different conditions, considering or not adjuvant and Ca^{2+} ions:

- 1) In-solution CaM;
- 2) In-solution CaM 10mM CaCl_2 ;
- 3) Adjuvant-adsorbed CaM;
- 4) Adjuvant-adsorbed CaM 10mM CaCl_2 .

Significant change in methionine oxidation rate slopes was observed in both cases of sample 2) and 4) respect to their no- Ca^{2+} reference samples 1), and 3).

This study demonstrated that HRF was able to detect induced conformational changes of CaM when exposed to Ca^{2+} ions in aluminum hydroxide-adsorbed samples.

Novel Aspect

HRF is a new application for interface-adsorbed protein structural characterization by MS. Implementation of this technique could provide an valuable analytical tool for adjuvanted vaccine efficacy and quality assessment.

References

1. R. A. Goldsby, T. Kindt and B. A., Immunology, 63-65 (2001).
2. C. Castelain, C. Genot, Biochim. Biophys. Acta, 1199, 59–64 (1994).
3. D. Kowalczyk, S. Slomkowski, F. W. Wang, J. Bioact. Compat. Polym., 9, 282-309 (1994).
4. M. C. L. Maste, W. Norde, A. J. W. G. Visser, J. Colloid Interface Sci., 196, 224-230 (1997).
5. F. C. Stevens, Can. J. Biochem. Cell Biol., 61, 906-910 (1983).
6. M. Chance, Wiley, (2008).

WP-24 / RAPID CHARACTERIZATION OF OVERPRODUCED PROTEINS USING NATIVE MASS SPECTROMETRY

Shay Vimer⁽¹⁾ - **Gili Ben-Nissan**⁽¹⁾ - **Hadas Cohen-Dvashi**⁽²⁾ - **Ely Morag**⁽¹⁾ - **Yoav Peleg**⁽³⁾ - **Tamar Unger**⁽³⁾ - **Meital Yona**⁽³⁾ - **Ron Diskin**⁽²⁾ - **Michal Sharon**⁽¹⁾

Weizmann Institute of Science, Department of Biomolecular Sciences, Rehovot (1) - Weizmann Institute of Science, Department of Structural Biology, Rehovot (2) - Weizmann Institute of Science, Structural Proteomics unit, Rehovot (3)

Keywords: native mass spectrometry, crude proteins samples, recombinant proteins, proteins analysis, ion mobility mass spectrometry

Introduction

Characterization of over-expressed proteins is usually coupled to costly and time-consuming purification steps. We previously presented a native MS approach for protein analysis directly from crude bacterial lysates without purification (1). Although prokaryotes expression systems are very efficient, they lack the machinery required for PTMs, highlighting the need for eukaryotic hosts. Here we present a method for characterization of proteins directly from crude samples from such systems.

Methods

As model systems we used different eukaryotic hosts that secrete recombinant proteins, the yeast *P. pastoris* expressing the cellulose binding module (CBM) from *C. thermocellum* and human serum albumin (HSA), and Sf9 insect cells expressing Transferrin receptor 1 (Tfr-1). Growth media from cells were collected and buffer exchanged into a MS compatible buffer. The different samples were analyzed using either an Orbitrap or a Synapt G2 HDMS QTOF-based platform (2).

Results

Initially we analyzed by native MS samples directly from *P. pastoris* and Sf9 growth media overcoming the need for purification procedure. Examination of the different growth conditions allowed us to monitor expression already after 24h of induction, from cultures, grown in 100 μ l medium in 96 well plates. We were able to monitor co-factor binding and validate that Ca²⁺ is associated with the CBM binding pocket using incubation of media with EGTA and calcium acetate. Furthermore, IM-MS spectra were measured for CBM directly from *P. pastoris* growth media. Using arrival time distribution, the CCS (Ω) was calculated and found to be highly correlated to theoretical CCS. Next, we examined the Tfr1 protein, secreted from Sf9 cells. Analysis of the crude growth medium generated highly resolved spectra, which revealed four different glycosylation states of the protein. We validated this by releasing the *N*-linked oligosaccharides by PNGase-F.

Conclusions

Here we show that direct native MS can be efficiently applied to different eukaryotic secretion systems. By examining the culture media of different expression systems we demonstrates that it is possible to rapidly determine the identity, solubility, oligomeric state, ligand binding, modification level and overall structure of secreted proteins. Moreover, we show that our method is not dependent on a specific MS platform but rather it is a general approach that can be applied on various instruments.

Novel Aspect

Here, we show that direct-MS can dramatically shorten the time gap between protein production and characterization.

References

1. Gan J., Ben-Nissan G., Arkind G., Tarnavsky M., Trudeau D., Noda Garcia Lianet., Tawfik D., Sharon M., *Analytical Chemistry*, 89, 4398-4404 (2017).
2. Hernández H., Robinson, C. V., *Nat. Protoc.* 2, 715–726 (2007)

WP-25 / INTEGRATIVE STRUCTURAL PROTEOMICS ANALYSIS OF THE 20S PROTEASOME COMPLEX

Rosa Viner⁽¹⁾ - *David M. Horn*⁽¹⁾ - *Eugen Damoc*⁽²⁾ - *Albert Konijnenberg*⁽³⁾

Thermo Fisher Scientific, San Jose (1) - Thermo Fisher Scientific, Bremen (2) - Thermo Fisher Scientific, research and development, Eindhoven (3)

Keywords: Native MS, crosslinking, top-down, pseudo MS3

Introduction

The structural determination of protein complexes plays an important role in the fundamental understanding their catalytic function. Multiple analytical methods, including hybrid approaches such as cryo-electron microscopy, cross-linking, and native mass spectrometry, are usually required for such type of analysis [1]. In this study we combined multiple mass spec based structural proteomics techniques to characterize the rabbit 20S proteasome complex.

Methods

Rabbit 20S proteasome complex was obtained from Boston Biochem. LC/MS bottom up, crosslinking and intact/top down analysis were performed using a Thermo Scientific™ UltiMate™ 3000 RSLC nano system and a Thermo Scientific™ Orbitrap™ Lumos™ mass spectrometer. Native MS experiments were performed using a modified Thermo Scientific™ Q Exactive™ Plus MS[2]. Data were analyzed with Thermo Scientific™ BioPharma Finder3.0, Proteome Discoverer™ 2.3, and ProSightPC™.

Results

The 20S proteasome complex contains 28 subunits arranged into four stacked rings: seven alpha non catalytic subunits and seven beta-subunits[3]. In the first series of experiments we performed bottom up, intact and top down analysis using databases of proteasome proteins from rabbit, human, and all species in UniProt. We identified 224 proteins and 17 different types of proteasome subunits and used the results to verify or correct the originally homology-based rabbit sequences. We then used this curated rabbit 20S proteasome database for all other experiments. For structural characterization, we performed native MS intact and top down experiments using a modified Q Exactive Plus mass spectrometer. The measured mass of the 28 subunits complex was consistent with expected value of 716.4 kDa and pseudo-MS3 experiments enabled identification of alpha 2 and beta 6 subunits. Crosslinking experiments using DSSO are combined with restraints from native and top-down MS to validate the structure obtained through electron microscopy.

Conclusions

Combining results from different MS methods, including novel native top-down -pseudo MS3 approach, we were able to perform comprehensive analysis of the rabbit 20S proteasome.

Novel Aspect

Complete workflow for structural characterization of protein complex

References

1. Greber, B. J. et al, Science, 348, 303-308, 2015
2. Van de Waterbeemd , M. et al, Nat. Methods, 14, 283-286, 2017
3. Bochtler, M et al, Ann. Rev. Biophys. Biomol. Struct. 28, 295-317, 1999

WP-26 / ASSEMBLY STUDIES OF ACETOBACTERIUM WOODII F1FO-ATP SYNTHASE AND ANALYZING CONFORMATIONAL CHANGES OF TMRABBY NESI

Khanh Vu Huu - *Stefan Brüchert* - *Robert Tampé* - *Volker Müller* - *Nina Morgner*

Wednesday

Institute of Physical and Theoretical Chemistry, Goethe-University, Frankfurt am Main, Germany

Keywords: Assembly of F1, LILBID-MS, ABC transport complex, ion mobility

Introduction

The F1FO ATP synthase of *Acetobacterium woodii* and the heterodimeric ABC transport complex TmrAB from *Thermus thermophilus* are responsible for ATP synthesis/hydrolysis and ATP-dependent membrane transport mechanism, respectively. We are interested in the assembly of subunits and investigate requirements for the different assembly steps. For TmrAB we apply ion mobility measurements to analyze conformational changes in dependency of different ligands.

Methods

In LILBID-MS a droplet generator generates droplets of an aqueous sample. These are transferred into vacuum and irradiated with an IR laser ($\lambda = 2.94 \mu\text{m}$) which leads to an explosive expansion of the droplets, ionization and finally a detection by a time-of-flight analyzer [1]. The Synapt G2-S HDMS offers ESI-mass spectrometry combined with ion mobility (IM) measurements.

Results

First we focused on the soluble part F1 of the F1FO ATP synthase of *Acetobacterium woodii*. The plasmids pKV003 (7036 bp), pKV004 (6928 bp) and pKV019 (6822 bp) (which include either the *atpA*- (encodes α), the *atpD*- (encodes β) or *atpG-atpC*(encodes $\gamma\epsilon$) genes were transformed in competent *E. coli* BL21gold (DE3) cells. The α -, β -subunit and the $\gamma\epsilon$ -complex were successfully purified by tandem affinity (Ni-chelating and Strep-Tactin) chromatography and can be characterized by LILBID-MS. Recently we are able to show the assembly efficiency of α -, β -subunits and also with the central stalk $\gamma\epsilon$ -complex of the ATP synthase in dependency of ATP and Mg^{2+} . Furthermore the steps of assembly were investigated by varying the laser energy.

A second project investigates three different ion mobility experiments we monitored conformational changes (inward/outward facing structure) of the heterodimeric ABC transport complex TmrAB related to ligand transport over the transmembrane by adding ATP/ADP at permissive conditions[2][3].

Conclusions

At soft LILBID conditions we could show that the assembly efficiency in the soluble part F1 can be improved in an ATP/ Mg^{2+} environment. Unfortunately the assembly is not yet efficient as we expected. We will vary the position for the purification tags to improve the assembly with the central stalk $\gamma\epsilon$.

With ion mobility experiments we are now able to monitor conformational changes of the ABC transport complex in dependence of substrate binding and addition of nucleotides.

Novel Aspect

Investigation of the assembly of a protein complex, via native MS studies. Analyzing conformational changes of membrane proteins with ion mobility.

References

1. Morgner N., Barth H. D., Brutschy B., Aust. J. Chem., 59, 109-114 (2006).
2. Bechara C. ; Nöll A. ; Morgner N. ; Degiacomi M. T. ; Tampé R., Robinson C. V., Nature Chemistry, 7, 255-262 (2015).
3. Zutz A., Hoffmann J., Hellmich U. A., Glaubitz C., Ludwig B., Brutschy B., Tampé R., J. Biol. Chem., 286, 7104-7115 (2010).

WP-27 / GAS-PHASE ELECTROPHORESIS (nES GEMMA) OF VIRUS-LIKE PARTICLES: RELATING PARTICLE SURFACE-DRY SIZE AND MOLECULAR WEIGHT

Victor U. Weiss⁽¹⁾ - ***Ronja Pogan***⁽²⁾ - ***Samuele Zoratto***⁽¹⁾ - ***Charlotte Uetrecht***⁽²⁾ - ***Günter Allmaier***⁽¹⁾

TU Wien, Inst. Chem. Technologies and Analytics, Vienna (1) - Heinrich-Pette-Institute (HPI), Hamburg (2)

Keywords: Gas-phase electrophoresis, Virus-like particle, VLP, nES GEMMA, native MS

Introduction

Gas-phase electrophoresis (on a nES GEMMA aka nES DMA) [1] separates surface dry particles according to their electrophoretic mobility diameter (EMD) which corresponds to the size of spherical nanoparticles. As shown for proteins and viruses [2,3], the EMD also relates to the molecular weight (MW) enabling to calculate an analyte MW based on its EMD. We now focus on virus-like particles (VLPs) and show the deviation of their EMD / MW correlation from other analyte cases.

Methods

VLP samples were desalted and transferred to ammonium acetate solution prior to native ESI MS and nES GEMMA analysis. Data for the following VLPs were combined / obtained: Human Rhinovirus Serotype 2 (HRV-A2) empty subviral particles [3], Hepatitis B Virus VLPs with T3 and T4 symmetry [4] and several norovirus VLP strains.

Results

VLPs of sufficient purity and concentration can be analyzed via native ESI MS and nES GEMMA. Comparing EMD values from nES GEMMA and MW values from native ESI MS experiments yields an EMD / MW correlation for VLPs. The particle MW is related to the EMD via a function to the power of approx. 3.5. Therefore, the EMD / MW correlation for VLPs was found to be significantly different from a similar correlation obtained for viruses (power of approx. 3.7) [3] and proteins (power of approx. 2.9) [5]. Hence, in order to be able to calculate reasonable MW values based on EMD data, the substance class (protein, intact virus, VLP) of the analyte in question needs to be known. If this requirement is met, nES GEMMA based calculation of a particle MW represents an interesting alternative at atmospheric pressure to the native MS approach.

Conclusions

We demonstrate that an EMD / MW correlation for VLPs can be obtained. Based on this data and EMD values obtained via nES GEMMA, MW values for VLPs of unknown MW can be calculated. In addition, we demonstrate the deviation of EMD / MW correlations for proteins and virus particles from the VLP case. A larger VLP dataset will further improve our EMD / MW correlation for VLPs.

Novel Aspect

We (i) present an EMD / MW correlation for VLPs and (ii) demonstrate the deviation of this VLP correlation from the protein and intact virus case.

References

1. Kaufman, S. L., Skogen, J. W., Dorman, F. D., Zarrin, F., Lewis, K. C., Anal. Chem. 68, 1895 (1996).
2. Bacher, G., Szymanski, W. W., Kaufman, S. L., Zollner, P., Blaas, D., Allmaier, G., J. Mass Spectrom. 36, 1038 (2001).
3. Weiss, V. U., Bereszczak, J. Z., Havlik, M., Kallinger, P., Gosler, I., Kumar, M., Blaas, D., Marchetti-Deschmann, M., Heck, A. J., Szymanski, W. W., Allmaier, G., Anal. Chem. 87, 8709 (2015).
4. Bereszczak, J. Z., Havlik, M., Weiss, V. U., Marchetti-Deschmann, M., van Duijn, E., Watts, N. R., Wingfield, P. T., Allmaier, G., Steven, A. C., Heck, A. J., Anal. Bioanal. Chem. 406, 1437 (2014).
5. Weiss, V. U., Golesne, M., Friedbacher, G., Alban, S., Szymanski, W. W., Marchetti-Deschmann, M., Allmaier, G., Electrophoresis, in press (2018)

WP-28 / NATIVE MASS SPECTROMETRY TO IDENTIFY STOICHIOMETRIES AND STABILITY OF CRISPR-CASCADE COMPLEXES

Sabine Wittig⁽¹⁾ - **Marius Rutkauskas**⁽²⁾ - **Ralf Seidel**⁽²⁾ - **Carla Schmidt**⁽¹⁾

HALOmem, Charles Tanford Protein Center, Institute for Biochemistry and Biotechnology, Martin Luther

University, Halle (1) - Molecular Biophysics group, University Leipzig, Leipzig (2)

Keywords: Native mass spectrometry, protein stoichiometry, CRISPR-Cas complexes

Introduction

The CRISPR-Cas system of bacteria or archaea is an adaptive immune system against foreign DNA. The type I system uses a large ribonucleoprotein surveillance complex called 'Cascade' in which the crRNA targets the foreign DNA by base-pairing and thus triggers DNA degradation. We used native mass spectrometry (MS) to determine the protein stoichiometries of CRISPR-Cascade complexes containing longer or shorter crRNAs than the wild-type complex.

Methods

Streptococcus thermophilus (St) CRISPR-Cascade was purified as described [1]. The purification buffer was exchanged against 200 mM ammonium acetate prior to native MS. Intact St-CRISPR-Cascade complexes were then analyzed on a Q-ToF Ultima mass spectrometer modified for transmission of high mass complexes [2] using gold-coated capillaries prepared in-house [3]. Mass spectra were processed and analyzed using MassLynx and MassSign software [4].

Results

St-CRISPR-cascade complex was purified and cascade subunits were identified by LC-MS/MS and database searching. All five subunits (*i.e.* Cse1, Cse2, Cas5, Cas6 and Cas7) were identified.

Next, we used native MS to define the stoichiometric of the wild-type CRISPR-cascade complex containing crRNA, Cse1, Cse2, Cas5, Cas6 and Cas7 subunits. The native MS spectrum revealed a Cse11(Cse2)2(Cas7)6-Cas51Cas61crRNA1 complex. The native MS spectrum also showed subcomplexes which lost Cse1 and Cse2 subunits as well as some copies of Cas7.

We further analyzed the impact of shorter or longer crRNA on the stoichiometry of the cascade complex. For this, we employed crRNAs which are 12 nucleotides longer or shorter. The native MS spectrum of the complex containing longer crRNA showed one abundant CRISPR-cascade complex with a stoichiometry of Cse11(Cse2)3(Cas7)8Cas51Cas61crRNA1. In contrast, the mass spectrum of complex containing shorter crRNA revealed several complexes; the most abundant complex was Cse11(Cse2)1(Cas7)4-Cas51Cas61crRNA1.

Conclusions

We conclude that an increase/decrease of 12 nucleotides in the crRNA results in addition or loss of two Cas7 and one Cse2 subunits. The CRISPR-cascade complexes containing longer crRNA appear to be more stable when compared with the shorter crRNA containing complexes and form mostly one abundant complex. Currently, we analyse Cascade complexes with variable crRNA length to gain insights into their stoichiometries and its impact on complex stability.

Novel Aspect

We studied the impact of crRNA length on protein stoichiometry and complex stability of CRISPR-Cascade complexes.

References

1. Sinkunas T, Gasiunas G., Waghmare S. P., Dickman M. J., Barrangou R., Horvath P., Siksnys V. *EMBO J.*, 32, 385-394 (2013).
2. Sobott F., Hernandez H., McCammon M.G., Tito M.A., Robinson C.V. *Anal. Chem.* 74, 1402-7 (2002).
3. Hernández H., Robinson C. V. *Nat. Protoc.* 2, 715-726 (2007).
4. Morgner N., Robinson C. V. *Anal. Chem.* 84, 2939-2948 (2012).

WP-29 / CHARACTERIZATION OF DEUTERATED ALPHA-CRYSTALLINE USING BENCH-TOP MALDI-TOFMS

Yuzo Yamazaki⁽¹⁾ - Ken Morishima⁽²⁾ - Nobuhiro Sato⁽²⁾ - Rintaro Inoue⁽²⁾ - Masaaki Sugiyama⁽²⁾

Shimadzu Corp., Global application development center, Kyoto (1) - Research Reactor Institute, Kyoto University, Osaka (2)

Keywords: MALDI-TOFMS, deuterated protein, crystalline, cross-linker

Introduction

Deuteration of protein is one of the important techniques in structural biology by using NMR or neutron scattering method. Incorporation of deuterated subunit protein into the hydrogenated one brings a new aspect of molecular dynamics in formation of multimer. As recent techniques for the deuteration have been more precise and more cost-effective, a quick and facile method to characterize deuterated proteins has been anticipated. In this report, we applied newly developed bench-top MALDI-TOFMS to analysis of deuterated β -crystalline to study substitution of hydrogen to deuterium in three different ways, molecular weight analysis, in-source decay, and tryptic peptide analysis. Furthermore, we will report multimer formation of the protein by using high mass detector and cross-linking chemistry.

Methods

Human recombinant hydrogenated/deuterated β -crystalline was prepared with overexpression in *E. Coli* with/without deuterated culture medium, as described[1]. Three samples, deuteration ratios of which were roughly estimated as 0%, 73%, and 98% (denoted as D0, D73, and D98, respectively), were subject to MALDI-TOFMS. Concentrations of them were 3 to 15 mg/mL. Bench-top MALDI-TOFMS (MALDI-8020, Shimadzu Corp., Japan) was applied to perform molecular weight measurement of whole proteins, in-source decay, and tryptic digest analysis in linear mode. Sinapinic acid, 1,5-DAN, and CHCA were used in those experiments. A cross-linker reagent was obtained from Covalx (Switzerland). All samples were reacted with the linker following an instruction provided from the company. After the reaction, desalted samples were mixed with sinapinic acid, then applied onto a stainless MALDI plate. MALDI-TOFMS (AXIMA-Performance, Shimadzu Corp., Japan) equipped with HM3 high mass detection system (Covalx, Switzerland) was applied to detection of cross-linked multimer in positive linear mode. The TOF analyser was calibrated using BSA and its multimer.

Results

Even though the newly developed bench-top MALDI-TOF MS has only linear mode, its performance was sufficient to estimate the deuteration ratio of the protein. Whole protein analysis of three samples was conducted with sinapinic acid, and mass shifts against D0% (hydrogenated) protein were estimated. While theoretical shift for 100% deuterium substituted protein was 1107.5, an observed shift of D98% was 1100, which indicates 99% deuteration ratio as almost expected. As well, an observed mass shift of D73% was 762.8, indicating 69% deuteration ratio, which was a little lower than expected. In-source decay using 1,5-DAN was conducted to examine an internal sequence of the proteins. Despite lower mass resolution caused by higher laser irradiation, a series of c-ions, observed up to 67 residues, were well-assigned to deuterated amino acids within the modest error, suggesting almost no significant bias of the deuterium substitution at any amino acid. Then, three proteins were digested with trypsin and isotopical distributions of each digested peptides were investigated. By simulation of isotopical distribution, deuteration ratios of the peptides were estimated as 99%, which was consistent with the value obtained with whole protein analysis. Additionally, hydrogenated β -crystalline was subject to cross-linking MS with high mass detector to observe multimer formation. As preliminary result, widely dispersed multimers were detected at the mass range from m/z 200K to 500K. Interestingly, whereas an intense monomer signal and a weak dimer signal were also observed, there was no intermediate formation of multimer, which is consistent with the proposed concept of "traveling subunits" that could be transiently liberated in subunit exchange[1].

Conclusions

Deuteration of protein was investigated with a newly developed bench-top MALDI-TOFMS. Molecular weight analysis, digested peptide analysis, and ISD enabled to estimate a deuteration ratio, which is a critical factor in study of a dynamics of protein multimer with small angle neutron scattering. The accessible

bench-top instrument could contribute to structural biology in prior use for the neutron scattering or NMR. Furthermore, the multimer of the crystalline was stabilized with cross-linking chemistry, and observed with high mass detector, through-put of which could be higher than native MS using ESI in analysis of multiple subjects.

Novel Aspect

A deuterated ratio of ^2B crystalline was estimated with a newly developed bench-top MALDI-TOFMS, and multimer of the protein was observed with cross-linking chemistry and high mass detector successfully.

References

1. Takata T, et al., Scientific Reports, 6, 29208 (2016)

WP-30 / DEVELOPMENT OF A CRYO-STAGE FOR LESA MASS SPECTROMETRY-TOWARDS TRULY NATIVE SURFACE SAMPLING OF PROTEINS

Bin Yan - Adam Taylor - Josephine Bunch

National Physical Laboratory, National Centre of Excellence in Mass Spectrometry Imaging, Teddington, UK

Keywords: Liquid Extraction Surface Analysis, Cryo sampling, Native mass spectrometry, Ion mobility spectrometry, Protein unfolding/refolding,

Introduction

Liquid extraction surface analysis (LESA) is a powerful method for measuring proteins from surfaces. In this work we present development and initial testing of a cryo platform for LESA MS. We explore the use of native sampling solutions for probing proteins directly from frozen surfaces.

Methods

An Advion Triversa Nanomate nano-electrospray source coupled with a Waters Synapt G2Si mass spectrometer was used for LESA MS. Modifications to the sample recess were made to incorporate water cooled copper bars and a thermoelectric peltier cooling system. Protein and LESA sampling solutions were prepared using either 100 mM ammonium acetate or 1:1 water/methanol with 0.1% formic acid. Experiments were conducted at room temperature and using the frozen stage.

Results

Preliminary results were obtained for ubiquitin. In ammonium acetate (20°C), native LESA sampling yielded ions with a narrow charge state distribution (5+ and 6+). However, H₂O/MeOH with 0.1% FA as the solvent resulted in a broad charge state distribution of ubiquitin ions (5+ to 13+) and slower average ion-mobility, indicating protein unfolding. For ubiquitin dissolved in H₂O/MeOH with 0.1% FA, denatured sampling also produced ions with a broad charge distribution. However, surprisingly, the mass spectrum acquired by native sampling of this denatured protein is dominated again by only 5+ and 6+ ions, which could be the evidence of protein re-folding. This was further supported by faster average mobility. The results using cryo sampling show very similar trend at room temperature, though the relative abundances of +5 and +6 ions differed for native LESA of denatured vs native proteins (1:2 for the former; 2:1 for the latter). The ion mobility was slightly slower than that of the original folded protein ion.

Conclusions

Protein unfolding or re-folding was observed to occur in LESA mass spectrometry study. The charge state distributions in the acquired mass spectra can not necessarily be used to judge whether the original structure of the studied protein is folded or unfolded. By employing cryo sampling method, the extent of protein unfolding or re-folding during LESA process can be largely minimized, therefore making truly native

surface sampling of protein promising.

Novel Aspect

A cryo-stage was developed for LESA-MS. With this protein re-folding and unfolding can be minimized during sampling so truly native surface analysis could be fulfilled.

WP-31 / NATIVE MASS SPECTROMETRY ANALYSIS OF NOROVIRUS GLYCAN INTERACTION

Hao Yan⁽¹⁾ - **Julia Lockhauserbäumer**⁽¹⁾ - **Katharina Bücher**⁽²⁾ - **Robert Creutzmacher**⁽³⁾ - **Alvaro Mallagaray**⁽³⁾ - **Thomas Peters**⁽³⁾ - **Laura Hartmann**⁽²⁾ - **Charlotte Uetrecht**^(1,4)

Leibniz Institute for Experimental Virology, Heinrich-Pette-Institut, Hamburg (1) - Institute of Organic and Macromolecular Chemistry, Heinrich-Heine-University Düsseldorf, Düsseldorf (2) - Institute of Chemistry, University of Lübeck, Lübeck (3) - European XFEL, Heinrich Pette Institute, Leibniz Institute for Experimental Virology, Hamburg (4)

Keywords: Native Mass Spectrometry, Noroviruses, Glycan, P domain, Binding.

Introduction

Human noroviruses (HuNoV) are non-enveloped RNA viruses belonging to the Caliciviridae family and the main cause of viral gastroenteritis. From the crystallographic structure, it is apparent that the norovirus capsid protein protruding (P) domain binds to glycans on the cell to mediate cell attachment. A detailed investigation of P domain-glycan interaction by native mass spectrometry (MS) is used to shed light on the mechanism behind viral cell entry.

Methods

Native MS is applied to investigate glycan binding to different HuNoV strains. Blood group antigen B tetrasaccharide was used to characterize glycan binding affinity and build titration curves. Other groups of carbohydrates such as gangliosides, human milk oligomers and globosides have been measured. We also tested multivalent glycan mimetics with different spacing and number(s) of glycan of interest for binding to P dimers to resemble glycan-binding in the natural situation. The number of binding events per type of natural saccharides and multivalent ligands will be compared and investigated between norovirus genotypes to elucidate binding profiles of glycans.

Results

The results elucidate the dissociation constants (KD) of natural glycans and glycan mimetics in norovirus cell attachment. Blood group antigen B revealed highest binding affinity among tested glycan candidates with a range of 100 μ M KD for the first binding event. A similar value was obtained by Saturation transfer difference nuclear magnetic resonance (STD-NMR). Another galactose containing natural glycan however showed similar binding affinity in native MS, while NMR data shows that single galactose cannot bind. The glycan mimetics binding affinity increases with increasing number of galactoses. The fucose containing mimetics showed good affinity but no avidity yet indicating that the spacing or rigidity requires further optimization. Another set of glycan mimetics based on these results are currently under investigation.

References

1. Mallagaray A, et al. (2015). *Angew. Chem.* 54, 12014–12019.
2. El Hawiet A, et al. (2012). *Anal. Chem.* 84, 50–58.
3. Cao S, et al. (2007). *J. Virol.* 81, 5949–5957.
4. Prasad BVV, et al. (1999). *Science* 286, 287–290.
5. Sun J, et al. (2006). *Anal. Chem.* 78, 3010–3018.

WP-32 / TOP-DOWN AND MIDDLE DOWN SEQUENCING OF IMMUNOGLOBULIN USING ELECTRON CAPTURE DISSOCIATION – TIME OF FLIGHT MASS SPECTROMETER COMBINED WITH ONLINE DISULFIDE BOND REDUCTION

Takashi Baba - Yves J. C. Le Blanc - Pavel Ryumin - Bill Loyd

Sciex, Research, Concord, Canada

Keywords: electron capture dissociation (ECD), antibodies, top down analysis, middle down analysis, online disulfide bond reduction

Introduction

Intact antibody sequencing is one of the goals in current mass spectrometry. Electron capture dissociation (ECD) provides unique features on top-down sequencing of proteins, which can be an ideal tool to analyze intact antibodies. We have developed a small and high throughput ECD device based on an RF ion trap [1]. This technology was applied to an antibody. We also implement an online disulfide bond reduction technique [2].

Methods

The ECD cell [1] was installed between Q1 and Q2 in a research grade quadrupole-TOF system. Typical electron beam irradiation time was 10 ms, and the electron beam intensity was tuned to obtain appropriate dissociation efficiency. The mass resolution of the TOF-MS is 35,000-47,000, which resolved isotope patterns of fragments up to $Z \sim 30^+$. A desalting LC column (Waters) was used for desalting, online reduction [2], and LC separation.

Results

Humanized monoclonal IgG (NIST-mAb) was obtained from NIST. This mAb was analyzed by the LC-ECD-TOF mass spectrometer. De novo sequencing on the intact ECD spectrum obtained by a single LC run indicated three sequences, and two of them were matched to N terminal partial sequences of the variable parts in a light chain and a heavy chain appeared in the human genome. The intact ECD spectrum was further analyzed in top-down manner using the suggested full sequences (the full sequence is provided by NIST), where the data covered the variable parts of the light chain and the heavy chain in the mAb. ECD at 3 eV did not cleave the disulfide-bonded rings in the protein. To obtain nearly complete sequence coverages, IdeS digestion (Genovis) was applied. For the online reduction of disulfide bonds (ref. 2), DTT was injected to the intact mAb and the IdeS digest trapped on the desalting column for 1 min.

Conclusions

By the reduction, we obtained sequence coverages of 84.7 % for the light chain, 78.3 % for the variable part of the heavy chain (Fd'), and 84.7% for the fixed part of the heavy chain (scFc). Further, ECD indicated the glycosylation site and its mass in scFc and CID informed the glycan composition.

Novel Aspect

High throughput antibody analysis using LC-ECD TOF mass spectrometer combined with online disulfide bond reduction technique

References

1. Takashi Baba et al. Anal. Chem. 2015, 87, 785–792
2. Yves Le Blanc et al. ASMSconference 2016 ThP-564

For information please contact: takashi.baba@sciex.com

WP-33 / CATION EXCHANGE CHROMATOGRAPHY IN DISPLACEMENT MODE FOR ONLINE ANALYSIS OF INTACT PROTEINS BY MASS SPECTROMETRY

Siti Nurul Hidayah - Dennis Kroesser - Benjamin Dreyer - Hartmut Schlüter

University Medical Center Hamburg Eppendorf, Institute for Clinical Chemistry and Laboratory Medicine, Mass Spectrometric Proteomics, Hamburg, Germany

Keywords: Displacement Chromatography, Top Down Proteomics, Intact Proteins

Introduction

Combination of mass spectrometry with chromatography is the key to successful mass spectrometric top down protein analysis [1]. Ion exchange chromatography in its classical form is incompatible with MS due to high salt concentrations in the eluent [2]. Displacement chromatography offers opportunity to perform cation-exchange chromatography (CEX) without using salt. In this study we tested if proteins eluting on-line from a CEX column into an ESI-MS can be detected.

Methods

An equimolar protein mixture (Lysozyme, Ribonuclease A, Cytochrome C and Myoglobin) was loaded (1.5 mg) onto a CEX column until it was saturated, recognized by the detection of the breakthrough of the proteins. Adsorbed proteins were eluted by pulsed injection of the displacer spermine. First 90 pulses contained 10 mM spermine, next 50 pulses 25 mM. The HPLC system was run with 0.1% formic acid as solvent and eluting proteins were detected by an Orbitrap mass spectrometer.

Results

Lysozyme and Ribonuclease A were eluted from the column by first 90 pulses of spermine (10 mM). Three different charge states from each protein were chosen and monitored in the corresponding spectra. Over the pulses, lysozyme and ribonuclease A eluted with constant concentration. A separation of these two proteins was not observed. After the injection of 140 pulses of spermine (25mM) Cytochrome C co-eluted with both previous proteins. Only Myoglobin stayed on the column and was not displaced by spermine. It eluted by washing the column with 1M ammonium acetate.

Conclusions

By application of displacement chromatography mode an online coupling of ion-exchange columns with MS detection for intact proteins is possible. However, the affinity of spermine obviously was not strong enough for displacing myoglobin. Thus, finding a suitable displacer with higher affinity for displacing proteins from CEX is still needed. Once found, it should be possible to even separate individual proteoforms of proteins.

Novel Aspect

Separation and detection of intact proteins by online coupling of ion-exchange columns with MS detection by application of displacement mode chromatography.

References

1. Muneeruddin K., Nazzaro M., Kaltashov I. A., Anal. Chem., 87, 10138-10145 (2015).
2. Vlasak J., Ionescu R., Curr. Pharm. Biotechnol., 9, 468-481 (2008).

WP-34 / HIGH AFFINITY EPITOPES OF APTAMER COMPLEXES OF THE MULTI-DOMAIN PROTEIN C-MET REVEALED BY PROTEOLYTIC AFFINITY MASS SPECTROMETRY

Loredana Lupu⁽¹⁾ - Francesca Rinaldi⁽²⁾ - Yannick Baschung⁽¹⁾ - Maxim V. Berezovski⁽³⁾ - Michael Przybylski⁽¹⁾

Steinbeis Center for Biopolymer Analysis and Biomedical Mass Spectrometry, University of Rostock,

Proteome Center, Rüsselsheim, Germany (1) - Laboratory of Pharmaceutical Analysis, University of Pavia, Department of Drug Sciences, Pavia, Italy (2) - John L. Holmes Mass Spectrometry Facility, University of Ottawa, Department of Chemistry and Biomolecular Sciences, Ottawa, Canada (3)

Keywords: Aptamer, Epitope, SPR, Maldi-MS, C-Met

Introduction

C-Met is a glycosylated receptor tyrosine kinase of the hepatocyte growth factor. Upon ligand binding C-Met transmits intercellular signals by its unique multi-substrate docking site. However, pathophysiological activation of the C-Met pathway leads to tumorigenesis, schizophrenia and cardiomyocytes death [1]. Here we report a study of the C-Met interaction with two DNA aptamers that bind the target protein (C-Met) with high affinity and specificity

Methods

The two C-Met aptamers used are CLN0003 and CLN0004 where produced using SELEX procedure [2]. For epitope identification, the aptamers where immobilized on sepharose beads. The binding epitopes were identified by proteolytic extraction and excizon-MS [3]. C-Met was digested with two enzymes. The unbound peptides where washed and the elution was done by slight acidification. Affinity determinations were carried out with a SPR biosensor.

Results

According to the SPR affinity data the aptamers are interacting with the protein with two different binding evets. The KDs where calculated with tree different methods to validate the results.

In the case of CLN0003 interacting with the protein we obtained the following KDs: $KD1 = 10^{-6}$ M and $KD2 = 10^{-8}$ M. For CLN0004 the KDs are as follows: $KD1 = 10^{-7}$ M and $KD2 = 10^{-8}$ M

Epitope analysis showed that the protein has indeed two different binding sites for the CLN003 aptamer. The two regions for the binding sites of the C-Met protein to the CLN0003 aptamer are (524-543) and (557-568). (SEELSGTWTQQICLPAYK and LTICGWDFGFR).

In the analysis of the epitope for the CLN0004 aptamer the results show that two different binding sites occur as well, one in the region of (251-257) with the corresponding sequence: HTRIIRF and one in the region of (381-388) with the sequence: NSSGCEAR

Conclusions

Because the C-Met protein has been linked to so many types of cancers it is important to determine the inhibiting mechanisms [4]. The work presented here can be used both for biomarker discovery, by using the CLN0003 aptamer or for delaying the cancer development. Furthermore, we present the discovery of the exact extracellular region necessary for inhibiting the interaction between the natural receptor and C-Met.

Novel Aspect

By proteolytic-MS we found for the first time the interaction sites of the C-Met protein with two different C-Met aptamers.

References

1. Organ SL, Tsao M-S, Ther. Adv. Med. Oncol., Nov; 3: S7–S19. (2011)
2. Iqbal, A., Labib, M., Muharemagic, D., Sattar, S., Dixon, B. R., Berezovski, M. V., PLoS ONE, 10-9 (2015)
3. Ștefănescu R, Lupu L, Manea M, Iacob RE, Przybylski M, J. Pept. Sci., 24-1 (2018)
4. Boltz A1, Piater B, Toleikis L, Guenther R, Kolmar H, Hock B., J. Biol. Chem., 286-24 (2011)

WP-35 / ANALYSIS FOR REFORMING PROCESSES FOR GASOLINE SURROGATE FUEL-AIR PRE-MIXTURES BY NON-THERMAL PLASMAS USING ION ATTACHMENT MASS SPECTROMETRY

Daiki Asakawa⁽¹⁾ - ***Eiichi Takahashi***⁽²⁾ - ***Naoaki Saito***⁽¹⁾

AIST, Research Institute for Measurement and Analytical Instrumentation, Tsukuba (1) - AIST, Research Institute for Energy Conservation, Tsukuba (2)

Keywords: Oxygen atom, Ozone, Ketone, Aldehyde, Ab initio calculation

Introduction

Non-thermal plasma induced reforming of fuel-air mixtures can improve the thermal efficiency of internal combustion engines by achieving its super-lean-burn operation. To address the general mechanisms of non-thermal plasma induced gasoline reforming in air, we used iso-octane, n-heptane, toluene, trimethyl pentene and methyl cyclohexane, which are components of gasoline surrogate fuel, as the models.

Methods

A mixture of air with liquid hydrocarbons was reformed by non-thermal plasma in a flow reactor, generated by a dielectric barrier discharge (DBD), and then directly analyzed using ion attachment mass spectrometry (IAMS). The reactor temperature was kept at 393 K and the equivalence ratio of the pre-mixture was set to 0.5. The simulation of the fuel reforming was performed by DFT calculation with Gaussian 16 program.

Results

The IAMS mass spectra of the alkane/air mixture with DBD irradiation showed intense signals corresponding to ketone and aldehyde. The DBD irradiation of the alkane/air mixture produced oxygen atom, which abstract a hydrogen atom from alkane. As a result, alkanes were converted to alkyl hydroperoxide through alkyl radicals. The ketone and aldehyde were generated by fragmentation of alkyl hydroperoxide. As in the case of oxygen atom generation, DBD produced ozone, which selectively react with alkene to produce molozonide. The trimethyl pentene undergoes the cleavage of double bond by DBD irradiation, leading to the ketones. The DBD irradiation to toluene/air mixture provided the signal corresponding to the benzaldehyde. A hydrogen atom in toluene was abstract by oxygen atom. The produced benzyl radical was immediately converted to benzyl hydroperoxide by the reaction with oxygen molecule and subsequent fragmentation would produce benzaldehyde. The proposed fuel reforming processes are supported by the ab initio DFT calculations.

Conclusions

DBD irradiation of the fuel/air mixture produced oxygen atom and ozone, which react with fuel molecules. Hydrocarbon was converted to ketone and aldehyde by DBD induced fuel reforming. The presence of products due to DBD plasma irradiation, such as ketone, aldehyde, and the intermediates *e.g.*, alkyl hydroperoxide and molozonide, would contribute to the combustion enhancement resulted in the thermal efficiency and stability of lean internal combustion engines.

Novel Aspect

DBD induced fuel reforming of hydrocarbon provide ketone and aldehyde, which would contribute to improve the performance of lean combustion in internal combustion engines.

Acknowledgements

This work was supported by the Council for Science, Technology and Innovation (CSTI), Cross-ministerial Strategic Innovation Promotion Program (SIP), and "Innovative Combustion Technology" (Funding agency: JST).

WP-36 / ELUSIVE MONOFUNCTIONAL ADDUCTS OF CISPLATIN WITH NATURAL AMINOACIDS REVEALED BY IRMPD SPECTROSCOPY AND PHOTOFRAGMENTATION KINETICS

Davide Corinti⁽¹⁾ - *Roberto Paciotti*⁽²⁾ - *Cecilia Coletti*⁽²⁾ - *Nazzareno Re*⁽²⁾ - *Susanna Piccirillo*⁽³⁾ - *Philippe Maitre*⁽⁴⁾ - *Bruno Bellina*⁽⁵⁾ - *Perdita Barran*⁽⁵⁾ - *Barbara Chiavarino*⁽¹⁾ - *Maria Elisa Crestoni*⁽¹⁾ - *Simonetta Fornarini*⁽¹⁾

Sapienza - University of Rome, Dipartimento di Chimica e Tecnologie del Farmaco, Roma (1) - Università G. D'Annunzio, Dipartimento di Farmacia, Chieti (2) - Università di Roma "Tor Vergata", Dipartimento di Scienze e Tecnologie Chimiche, Roma (3) - Université Paris-Sud, Laboratoire de Chimie Physique, Orsay (4) - University of Manchester, Manchester Institute of Biotechnology, Manchester (5)

Keywords: IRMPD kinetics, isomerism, conformer discrimination, cisplatin, aminoacids

Introduction

Relevant intermediate species in the reaction path of cisplatin with aminoacids can be difficult to characterize due to the high reactivity of square planar complexes of Pt^{II} with nucleophiles. However, ESI-MS proved to be valuable in allowing to detect for the first time the primary complexes generated from direct substitution of one of the labile chloro ligands of cisplatin with aminoacids like His and Met, modeling cisplatin interaction with peptides and proteins.

Methods

The complexes of interest were brought to the gas-phase using ESI and mass-analyzed using a Paul ion trap instrument (Bruker Esquire 6000) or a FT-ICR mass spectrometer (Bruker APEX-Qe and 7T Solarix). Vibrational features of the mass-selected ions were assayed using IRMPD spectroscopy in both XH (X = C, N, O) stretching and fingerprint regions using an OPO/OPA benchtop laser system and the free electron laser (FEL) at the Centre Laser Infrarouge d'Orsay, respectively.

Results

The primary complexes of cisplatin with histidine and methionine, cis-[PtCl(NH₃)₂(L)]⁺ (L = His, Met), have been characterized using high-resolution MS and IRMPD spectroscopy. The IR spectrum of bare ions can be directly compared with calculated spectra, thus allowing the attribution of the assayed species to a particular conformer or a set of them. In the case of the cisplatin-histidine adduct, a complex mixture of isomers and conformers was hypothesized and finally resolved using IRMPD kinetics on conformer-specific vibrational modes. [1] Regarding the adduct with methionine, the comparison of the IRMPD spectrum with theoretical ones confirmed the preferred interaction of Pt with the thioether functionality of the aminoacid. However, when the complex of a structural isomer of the drug, transplatin, with Met was assayed, differences in the photofragmentation kinetics on a same, selected vibrational mode permitted to unveil the presence of a mixture of isomers presenting platinum coordination at either the amino or the thioether group. [2]

Conclusions

We show the ability of mass spectrometry to extract and isolate from solution species difficult to characterize using condensed-phase based techniques. Also, the coupling with IRMPD permitted to obtain spectroscopic information on the mass-isolated ions, proving the presence of complex isomeric mixtures. In this regard, the use of fixed-wavelength IRMPD kinetics permitted to discriminate and quantify isomers and conformers forming the gas-phase populations.

Novel Aspect

Monofunctional adducts of cisplatin with His and Met have been characterized for the first time. IRMPD kinetics were used to discriminate and quantify isomers and conformers.

References

1. Corinti D., De Petris A., Coletti C., Re N., Chiavarino B., Crestoni M. E., Fornarini S., ChemPhysChem, 16, 318 (2017).
2. Paciotti R., Corinti D., De Petris A., Ciavardini A., Piccirillo S., Coletti C., Re N., Maitre P., Bellina B., Barran P., Chiavarino B., Crestoni M. E., Fornarini S., Phys. Chem. Chem. Phys., 19, 26697(2017).

WP-37 / CHARACTERIZATION OF A NITROGEN-OPERATED DIELECTRIC BARRIER DISCHARGE IONIZATION SOURCE

Luzia Gyr⁽¹⁾ - Felix David Klute⁽²⁾ - Joachim Franzke⁽²⁾ - Renato Zenobi⁽¹⁾

ETH Zurich, ETH Zurich, Zurich (1) - ISAS Dortmund, ISAS Dortmund, Dortmund (2)

Keywords: ambient ionization techniques, nitrogen dielectric barrier discharge, fundamental, reactive species

Introduction

In the last decade, ambient ionization coupled to mass spectrometry generated enormous attention due to the fast and sensitive analysis with nearly no sample preparation. [1] However, the development of new sources dominated over the understanding of the ionization mechanism of already existing methods. In this work, the active capillary plasma ionization source was characterized to get a better understanding of the behavior of the reactive species.

Methods

The active capillary plasma ionization source contains a stainless steel capillary (ground electrode) which is separated by a glass capillary from a copper ring (sine-modulated high voltage). [2] Nitrogen or air was used as discharge gas. The plasma source is directly connected to the mass spectrometer. All optical measurements were performed AT the cross section of the plasma source using either UV-VIS spectrometer in the range of 200-850 nm or a ns-gated iCCD camera.

Results

The active capillary plasma ionization source based on nitrogen dielectric barrier discharge was characterized by optical emission spectroscopy. The two most abundant signals were the transitions of NO (γ system) and N_2 (2nd positive system) and in low abundance the transition of N_2^+ (1st negative system) and O_2^+ (Schuhmann band) was obtained. It is known that N_2^+ is the reactive species for producing $(H_2O)H_3O^+$ or SH^+ , which are mostly responsible for the protonation of the analytes (MH^+). However, the lower abundance of N_2^+ was obtained in nitrogen plasma than in helium, which is due to the different ionization pathway. Nevertheless, the sensitivity of the ionization was similar for both discharges, since different geometries to the MS and sample introduction into the plasma was used. The nitrogen plasma ignites as irregular filaments between the inner electrode and the glass, which was confirmed by using an iCCD camera. Analytes can therefore travel through the plasma-free regions inside the capillary without fragmentation.

Conclusions

The properties of the filaments were studied by monitoring the transitions of NO (γ system) and N_2 (2nd positive system) when applying different plasma parameters (voltage, frequency, humidity, etc.). Characterization of the plasma not only helps to improve the understanding of the ionization mechanism but also allows to optimize the performance of the plasma source when coupled with the mass spectrometer.

Novel Aspect

This is the first report on spectroscopic characterization of the active capillary plasma ionization source

which is used as a soft ionization source coupled to MS.

References

1. X. Ding, Y. Duan, *Mass Spectrom. Rev.*, 34, 449–473 (2015).
2. L. Gyr, J. -C. Wolf, J. Franzke and R. Zenobi, *Anal. Chem.*, 90, 2725-2731 (2018).

WP-38 / A TARGETED METABOLOMICS APPROACH FOR DIAGNOSIS OF INBORN ERRORS OF METABOLISM

Minnie Jacob^(1,2) - ***Abeer Malkawi***^(1,3) - ***Noor Albast***⁽⁴⁾ - ***Salam Al Bougha***⁽⁴⁾ - ***Andreas Lopata***⁽²⁾ - ***Majed Dasouki***⁽¹⁾ - ***Anas Abdelrahman***^(1,4,5)

Department of Genetics, King Faisal Specialist Hospital and Research Center, Riyadh, Saudi Arabia (1) - Department of Molecular & Cell Biology, Australian Institute of Tropical Health and Medicine, James Cook University, Townsville, Australia (2) - Department of Clinical Pharmacy, Jordan University of Science and Technology (JUST), Irbid, Jordan (3) - College of Medicine, Al Faisal University, Riyadh, Saudi Arabia (4) - Department of Chemistry, Memorial University of Newfoundland, St. John's, NL (5)

Introduction

Metabolome, the end product of the genome, can be studied through identification and quantification of small molecules. The global metabolome influences the individual phenotype through clinical and environmental interventions. (1, 2) Metabolomics has become an integral part of clinical research and allowed for another dimension of better understanding of disease pathophysiology. The clinical biochemistry laboratory routine workload (>95%) is based on small molecular identification, which can potentially be analyzed and discovered through metabolomics. (3, 4) However, multiple challenges in clinical metabolomics impact the entire workflow and data quality, thus the biological interpretation needs to be standardized for a reproducible outcome. The existing Mass spectrometry-based Newborn Screening assay has false discovery rate (FDR) that can be minimized once linearly integrated with other markers.

Methods

Herein, we introduce the establishment of a comprehensive targeted metabolomics method for a panel of 225 clinically relevant metabolites using Liquid chromatography-tandem mass spectrometry (LC-MS/MS) standardized for clinical research. The sensitivity, reproducibility and molecular stability of each targeted metabolite (amino acids, organic acids, acylcarnitines, sugars, bile acids, neurotransmitters, polyamines, and hormones) were assessed under multiple experimental conditions. The metabolic tissue distribution was determined in various rat organs. Furthermore, the method was validated in dry blood spot (DBS), the standard sample type for newborn screening, collected from patients known to have various inborn errors of metabolism (IEMs).

Results

The targeted LC-MSMS-based metabolomics method was developed for quantitation purpose, and the validity of this method was studied by following the United States Food and Drug Administration (USFDA) guideline. (5) This approach appears to be sensitive (LLOQ for 2-Ketobutyric acid is 1 nM) and robust (CV% ranges from 80-120%). The metabolic chemical stability and extraction recovery were evaluated in this approach using several biological matrices to control the method robustness in the routine analysis. The sensitivity of the method demonstrated differential and unique expression profiles in various rat tissues, where the overall metabolic expression varied, but was highest in the kidney tissue (35%) compared to the liver (25%), brain (17%), heart (15%) and muscle tissue (8%) for the same amount of tissue. Several optimizations and studies have been done on DBS samples such as relative metabolic expression compared to whole blood and serum, the normal cutoff for each expressed metabolite in this panel, and the differential expression compared to patient samples. 56 patient samples with known diagnosis were analyzed using this platform for clinical validation. This method was capable to produce the same diagnosis

as in the routine NBS platform using the same key markers, for instance, it gives the Phenylalanine/tyrosine ratio high in PKU patients compared to the others. Additionally, several other metabolites were found to be differentially expressed in these patients compared to control samples such as betaine, succinic acid, niacinamide, fructose 1, 6-bisphosphate. Several sugars (*i.e.* arabinose, ribulose) were overexpressed mainly in the disorders of mitochondrial metabolism (VLCAD and GA2 deficiency) compared to controls. In MMA (n=6) and PA (n=2) patients displayed distinct metabolomics profiles with minimal overlap including hyperglycinemia and related elevated glutamine levels. Interestingly, spermine appears to be markedly reduced in MMA compared to PA suggesting a preferential toxicity to the synthetic polyamine pathway in such patients.

Conclusions

There is a strong clinical need to develop and validate quantitative methods capable of measuring the largest number of metabolites subject to wide and rapid dynamic changes present in a biological system. In the present study, we demonstrated the effect of sample matrix on the expression of profiles of a large set of clinically oriented metabolites which may help to guide the choice of the optimal sample type and matrix in which it is collected. We have also shown the capability of this assay to detect common inborn errors of metabolism. In addition to its role as a second tier assay, this targeted panel is also potentially useful as a screening tool for a large group of inborn errors of metabolism

Novel Aspect

This panel has the ability to identify patients with IEMs and discover novel markers that minimize the assay FDR.

References

1. Dettmer K, Aronov PA, Hammock BD. Mass spectrometry Reviews. 26, 51-78 (2007)
2. Jacob M, Lopata AL, Dasouki M, Abdel Rahman AM. Mass Spectrometry Reviews. 1-18 (2017)
3. Wang Q, Sun T, Cao Y, Gao P, Dong J, Fang Y, et al. OncoTargets and Therapy 9, 1389-1398 (2016)
4. Chace DH, Kalas TA, Naylor EW. Clinical Chemistry 49, 1797-1817 (2003)
5. Garofolo F, Michon J, Leclaire V, Booth B, Lowes S, Viswanathan CT, et al. Bioanalysis 4:231-236 (2012)

WP-39 / STEROID BISCONJUGATES: UNEARTHING HIDDEN TREASURES OF STEROID METABOLISM

Malcolm McLeod - Christopher Fitzgerald

Research School of Chemistry, Australian National University, Canberra

Keywords: steroid bisconjugate, molecular di-anion, stable isotope label, MS fragmentation, ion loss,

Introduction

Doubly or bisconjugated steroid metabolites have been known for decades and typically studied by indirect GC-MS methods. Recently, the direct detection of steroid bis(sulfates) by LC-MS has been applied to doping in sport [1] and the prenatal diagnosis of inborn errors of steroid biosynthesis [2]. Access to labelled compounds for MS would increase our understanding of these fascinating compounds.

Methods

General methods of chemical synthesis were developed to prepare stable isotope ^{18}O -labelled sulfate esters and ^{13}C -labelled glucuronides, including selectively labelled steroid bis(sulfate) and bisglucuronide compounds for use as MS probes or internal standards. The ionisation and fragmentation of these reference materials have been studied in detail on a triple quadrupole mass spectrometer.

Results

Stable isotope ^{18}O -labelled sulfate esters and ^{13}C -labelled glucuronides were prepared from commercially available $\{^{18}\text{O}_4\}$ -sulfuric acid (96 atom %) and $\{^{13}\text{C}_6\}$ -D-glucose (>99 atom %). High levels of incorporation were observed with the conjugates suitable for use as MS probes or stable isotope labelled internal standards. Selectively labelled steroid bis(sulfate) and bisglucuronide compounds were used to study MS fragmentation pathways. Steroid bis(sulfate) compounds preferentially ionized as the doubly deprotonated molecule ($[\text{M}-2\text{H}]^{2-}$) and this precursor typically underwent ion loss of the HSO_4^- ion (m/z 97) to afford fragment ions of higher m/z ($[\text{M}-2\text{H}-\text{HSO}_4]^-$). Other fragmentation processes including ion loss of $^{\bullet}\text{SO}_3^-$ (m/z 80) and HSO_3^- (m/z 81), neutral loss of SO_3 (40 Da) and ion molecule "roaming" were observed and correlated with analyte structure.

Conclusions

Access to stable isotope labeled sulfate and glucuronide conjugates has unveiled the rich MS fragmentation behaviors of steroid bisconjugates. This knowledge paves the way for the development of new LC-MS approaches for the direct and untargeted detection and identification of these metabolites.

Novel Aspect

Synthesis of stable isotope labelled sulfate and glucuronide conjugates, selectively labelled bisconjugates and their use in MS fragmentation studies.

References

1. McLeod M. D., Waller C. C., Esquivel A., Balcells G., Ventura R., Segura J., Pozo Ó. J., *Anal. Chem.*, 89, 1602–1609 (2017).
2. Pozo O. J., Marcos J., Khymenets O., Pranata A., Fitzgerald C., McLeod M. D., Shackleton C., *J. Mol. Endocrinol.*, DOI: 10.1530/JME-17-0286 (2018).

WP-40 / INVESTIGATING THE REACTION MECHANISM OF THE OXIDATION OF HEPTAFULVENES INTO TROPONES BY A MICROREACTOR COUPLED ONLINE TO MASS SPECTROMETRY AND DENSITY FUNCTIONAL THEORY CALCULATIONS

Sofia M. E. Nilsson⁽¹⁾ - *Henning Henschel*⁽²⁾ - *Gianmario Scotti*⁽¹⁾ - *Markus Haapala*⁽¹⁾ - *Alexandros Kiriazis*⁽¹⁾ - *Gustav Boije af Gennäs*⁽¹⁾ - *Tapio Kotiaho*^(1,3) - *Jari Yli-Kauhaluoma*⁽¹⁾

Drug Research Program, Division of Pharmaceutical Chemistry and Technology, Faculty of Pharmacy, University of Helsinki, P.O. Box 56 (Viikinkaari 5 E), FI-00014, Finland (1) - Research Unit of Medical Imaging, Physics and Technology, University of Oulu, P.O. Box 8000, FI-90014, Finland (2) - Department of Chemistry, Faculty of Science, University of Helsinki, P.O. Box 55, FI-00014, Finland (3)

Keywords: Reaction monitoring, density functional theory, microreactor, online, reaction mechanism

Introduction

Understanding reaction mechanisms is crucial in studies of chemical reactions. When knowing the mechanism, reaction conditions can be designed to maximize the yield of desired product(s). Microreactors coupled online to mass spectrometry are ideal for experimental reaction mechanism studies [1],[2], while density functional theory calculations are well-established for theoretical studies.

Methods

A 3D printed polypropylene microreactor with an integrated ESI needle [3] was interfaced with an ion trap mass spectrometer for online analysis of the oxidation of heptafulvenes to tropones. Density functional theory calculations were done with the software Gaussian 09, using the M06-2X functional and a 6-311++G(d,p) basis set to identify the most favored reaction pathway.

Results

Online mass spectrometry allowed observation of the heptafulvene starting material, two key intermediates of the oxidation reaction and the tropone product as protonated molecules in positive mode ESI. MS_n analysis of these species, as well as how their intensity change in relation to each other depending on the reaction time, support the proposed identifications. Based on the mass spectrometric results six competing reaction pathways were studied with density functional theory. The lowest energy reaction pathway, which thus is the most likely one, proceeds via initial formation of an epoxide, which is opened upon the addition of a second equivalent of the oxidizing species meta-chloroperoxy benzoic acid. The adduct then undergoes a Criegee-like rearrangement to produce a positively charged hemiketal, which upon deprotonation dissociates into acetone and the tropone.

Conclusions

The most likely reaction mechanism for oxidizing heptafulvenes to tropones has been identified by experimental and theoretical investigations. These results are a good basis for optimizing reaction conditions to improve yields of the promising anti-cancer compound tropone [4]. Furthermore, these results highlight the usefulness of combining experiments and theory for studying reaction mechanisms.

Novel Aspect

The synergy of combining results from a microreactor coupled to mass spectrometry and computational chemistry to find a plausible reaction mechanism.

References

1. Santos, L. S. Reactive Intermediates: MS Investigations in Solution, Wiley-VCH Verlag GmbH & Co. KGaA (2009).
2. Sikanen, T. et al., Mass Spectrometry Reviews, 29, 351–391 (2010).
3. Scotti, G. et al., Reaction Chemistry and Engineering, 2, 299–303 (2017).
4. Kiriazis, A. et al., PLoS One, 8, e55409 (2013).

WP-41 / HOW TO MODIFY THE ACIDITY OF CHARGED DROPLETS?

Mari Ojakivi⁽¹⁾ - ***Jaanus Liigand***⁽¹⁾ - ***Anneli Kruve***⁽²⁾

The University of Tartu, Institute of Chemistry, Tartu (1) - Freie Universität Berlin, Institut für Chemie und Biochemie, Berlin (2)

Keywords: aerosols, electrospray, ionization efficiency, protonation, surface chemistry

Introduction

The concept of acidity in confined spaces has been poorly understood, especially in media violating electroneutrality. Here [1] we present a study on the acidity of charged droplets to fill the gap in knowledge. We studied the effect of pH and different additives on the acidity of charged water and methanol droplets via protonation of nitrogen bases in electrospray. As a result, the work gives insight into how to modify the acidity of charged droplets.

Methods

The protonation of 10 nitrogen bases was studied using ESI/MS. The extent of protonation (which is related to the acidity of charged droplets) was calculated using ionization efficiency values. To alter the acidity of charged droplets produced in ESI, 28 solvent compositions were used in water (pH H₂O=1.42–10.29) and 14 in methanol (pH MeOH=3.63–12.24). The solvent composition was varied using 10 additives in different concentrations and combinations.

Results

Generally, we observed that the ionization efficiencies of most analytes are pH dependent in specific pH range in both water and methanol. However, wrong-way-round ionization occurs in both solutions. Also, the extent of protonation decreases sharply in the pH H₂O range from 2 to 4 for all compounds, even though pK_a H₂O values of these compounds range from 2.7 to 9.0. Close investigation of the pH range of 2 to 4 revealed that additive type affects the acidity of the droplets more than pH in aqueous solutions, enhancing the extent of protonation of the compounds when only acidic additives are present compared to droplets of the same pH containing also ammonium salts. Thus, we propose that in case of acidic additives, the excess charge on the droplet surface is a result of protonated acid and water molecules and, therefore, enhances ionization efficiency of compounds significantly. The effect of additives is less influential in case of methanol droplets, as the excess charge is situated in the interior of the droplets.

Conclusions

The results show that the concept of acidity in solution phase cannot be transferred to charged droplets. We demonstrated that the acidity of charged water droplets depends on the additive; which means that the extent of protonation can effectively be modified by changing additive type even if the pH is kept constant. The opposite is true for charged methanol droplets where the extent of protonation can be adjusted with the conventional approach of changing the pH.

Novel Aspect

While the protonation in charged methanol droplets depends on the solution-phase acidity, the protonation in charged water droplets can easily be modified with additives.

References

1. Ojakivi M., Liigand J., Krüge A., *Chemistry Select*, 3, 335–338 (2018)

WP-42 / STUDIES OF A L-PROLINE CATALYZED INVERSE ELECTRON DEMAND DIELS-ALDER REACTION BY ESI-MS

Anne Schnell - J. Alexander Willms - Marianne Engeser

University of Bonn, Research Group Engeser, Bonn

Keywords: Charge-tagging; electrospray ionization mass spectrometry; L-proline; organocatalysis; reaction mechanism

Introduction

The synthetic scope of the L-proline catalyzed inverse electron demand Diels-Alder reaction between tetrazines and ketones has been published by Xie et al. [1] in 2008 alongside a postulated catalytic cycle. The postulated catalytic cycle includes three intermediates. First an enamine I derived from L-Proline and a ketone, which undergoes the Diels Alder reaction with the tetrazine forming II which loses nitrogen in a retro Diels-Alder reaction forming III. [1]

Methods

Samples of the reaction solution were analyzed by ESI-MS and CID experiments. A charge-tagged substrate and a catalyst were utilized in different experiment setups as to enhance the ESI response factors of the respective intermediates.

Results

While studying the reaction without implementing a charge tag into any of the participating species the third intermediate III could be detected and kinetic studies could be obtained. The charge tagged tetrazine

was synthesized in hopes of detecting the second intermediate II as well, but to no avail; however again kinetic studies could be obtained. By using the charge tagged L-proline derivative from our group[2] as a catalyst all three intermediates could be detected and characterized by CID experiments. Within the CID experiments it was possible to mimic the reaction steps of the catalytic cycle.

Conclusions

It was possible to thoroughly verify the catalytic cycle for the L-proline catalyzed inverse electron demand Diels-Alder reaction between acetone and 3,6-di-2-pyridyl-1,2,4,5-tetrazine.

Novel Aspect

The superior effectiveness of utilizing a charge-tagged catalyst in ESI-MS studies of catalytic reactions could be demonstrated.

References

1. H. Xie, L. Zu, H. R. Oueis, H. Li, J. Wang, W. Wang, *Org. Lett.*, 10, 1923 (2008).
2. J. A. Willms, R. Beel, M. L. Schmidt, C. Mundt, M. Engeser, *Beilstein J. Org. Chem.*, 10, 2027 (2014).

WP-43 / FRAGMENTATION OF ARSENIC-CONTAINING CHEMICAL WARFARE AGENTS IN GC-EI/MS AND LC-ESI/MS

Martin Söderström - Paula Vanninen

VERIFIN, University of Helsinki, Helsinki

Keywords: organoarsenic chemicals, fragmentation pathways, rearrangements

Introduction

Arsenic-containing chemical warfare agents were used in World War I and World War II. Large part of the Axis stockpiles were destroyed after the World War II by burying in soil and dumping in sea. The environmental effect of these actions are still seen around the world. Mass spectrometry is widely used to analyze these chemicals and their degradation products. Arsenic atom induces a large number of rearrangement reactions.

Methods

The fragmentation pathways of selected organoarsenic chemicals were studied using gas chromatography-electron ionization/high resolution mass spectrometry (GC-EI/HRMS) and liquid chromatography-electrospray ionization/high resolution mass spectrometry (LC-ESI/HRMS). The latter instrument can produce also MSn/HRMS data using quadrupole, ion trap and Orbitrap analyzers. The spectra were carefully analyzed to identify fragmentation pathways.

Results

Some characteristic fragmentation pathways were identified based on the mass spectra. Arsenic atom induce rearrangement reactions, which lead formation of large number of cyclic fragmentation products. The chemicals discussed in this paper will include intact warfare agents, which are typically chlorine-containing chemicals, their degradation and oxidation products as well as some alkylthiol derivatives required for GC-based analyses.

For most of the chemicals there are several fragments, where neutral arsenic species are lost from the ionized molecule. One of the most exotic rearrangements is the formation of a phenyl cation from tris(2-chlorovinyl)arsine ("Lewisite I"). It is worth noting that in LC-ESI/MS, diphenylarsenic chemicals are producing highly characteristic radical cations, *e.g.* $C_{12}H_{10}^{+•}$, $C_{12}H_8^{+•}$ and $C_6H_5As^{+•}$.

The current paper will present the suggested fragmentation pathways and highlight common ions formed

from phenyl and diphenylarsenic chemicals.

Conclusions

Organoarsenic chemicals are an interesting group of chemicals, which fragment in unconventional ways to produce characteristic fragments. These fragment ions can be identified and then used for search of novel chemicals possibly forming in the environment. Several novel chemicals have been tentatively identified using these fragments in sediment samples taken in CHEMSEA, MODUM and DAIMON projects from dump sites in the Baltic Sea.

Novel Aspect

Novel chemicals have been detected in environmental samples. Their spectra have been analyzed using data obtained from other phenyl and diphenyl arsenic chemicals.

WP-44 / IDENTIFICATION AND INTERCONVERSION OF ISOMERIC 4,5-FUNCTIONALIZED 1,2,3-THIADIAZOLES AND 1,2,3-TRIAZOLES IN CONDITIONS OF ELECTROSPRAY IONIZATION

Marina Zimens⁽¹⁾ - *Dmitrii Mazur*⁽¹⁾ - *Vasiliy Bakulev*⁽²⁾ - *Albert Lebedev*⁽¹⁾

Lomonosov Moscow State University, Organic Chemistry Department, Moscow, Russian Federation (1) - Ural Federal University, Institute of Chemical Engineering, Ekaterinburg, Russian Federation (2)

Keywords: HRMS, electrospray ionization, 1,2,3-thiadiazoles, 1,2,3-triazoles, isomer identification, Orbitrap

Introduction

Investigation of 1,2,3-triazoles and 1,2,3-thiadiazoles has a great biological potential. However the problem of isomerisation of 4,5-functionalized 1,2,3-thiadiazoles and 1,2,3-triazoles in the presence of an acid or a base at room temperature makes it difficult to identify each isomer. The present study deals with ESI-MS application in positive and negative modes to identify each isomer by CID and HCD mass spectra and reproduce the rearrangement processes.

Methods

All experiments were performed with Orbitrap Elite mass-spectrometer (Thermo Fisher Scientific, USA) with electrospray ionization source. The acetonitrile solution of each compound was introduced through syringe pump directly into the ion source at 5 $\mu\text{l}/\text{min}$. High resolution mass spectra were acquired both in positive and negative modes. MS/MS experiments were carried out using both CID and HCD fragmentation techniques at different energies.

Results

The monomolecular transformations under (\pm)ESI-MS conditions of 4,5-functionalized 1,2,3-thiadiazoles and 1,2,3-triazoles in the gas phase were investigated varying the amino- and sulfanyl moieties. Fragmentation pattern of each isomer was carefully studied by scanning CID and HCD energies in positive and negative modes. It was demonstrated that basic fragmentation process in case of protonated 1,2,3-thiadiazoles involves the loss of N_2 molecule. All other fragment ions in MS2 spectra originate from the primary $[\text{MH}-\text{N}_2]^+$ ion. Unlike 1,2,3-thiadiazoles, MH^+ ions of 1,2,3-triazoles demonstrate a richer fragmentation pathway without losing N_2 molecule. The main fragmentation process of MH^+ ion involves the loss of sulfynylaryl/sulfynylalkyl or sulfonylaryl/sulfonylalkyl radical. Appearance of $[\text{MH}-\text{N}_2]^+$ ion in the spectra of 1,2,3-triazoles was proved to form due to the gas phase rearrangement during the ESI-MS analysis. The connection of all precursor and product ions was established in MS3 experiments.

Conclusions

Analysis of (\pm)ESI-MS/MS mass spectra of the isomeric 4,5-functionalized 1,2,3-triazoles and 1,2,3-thiadiazoles with HRMS allowed identifying each isomer. Moreover, only 1,2,3-triazoles undergo

transformation into the corresponding 1,2,3-thiadiazoles under (+)ESI-MS conditions, thus mimicking the reaction in the condensed phase. Alkyl substituent next to the sulfone group in both isomers leads to new fragmentation pathways.

Novel Aspect

1,2,3-Triazoles and 1,2,3-thiadiazoles have been differentiated using (\pm)ESI-MS/MS method. 1,2,3-Triazoles were shown to rearrange during the (+)ESI-MS analysis.

WPS-S03 - Pharmaceuticals

Chairs: Christine Eckers, Gianluca Bartolucci

WP-45 / ULTRASOUND ASSISTED EXTRACTION PARAMETERS OPTIMIZATION BY EXPERIMENTAL DESIGN FOR STEROIDAL SAPONINS FROM PLANT MATERIAL AND CELL CULTURES OF DIOSCOREA DELTOIDEA

Boris Sarvin⁽¹⁾ - **Timur Baygildiev**⁽¹⁾ - **Dmitry Kochkin**⁽²⁾ - **Andrey Stavrianidi**⁽¹⁾ - **Igor Rodin**⁽¹⁾

Lomonosov Moscow State University, Department of Chemistry, Moscow (1) - Lomonosov Moscow State University, Faculty of Biology, Moscow (2)

Keywords: ultrasound assisted extraction, steroidal saponins, Dioscoreadeltoidea plant cell culture, Latin square experimental design, HPLC-MS

Introduction

Many of steroidal saponins have anticancer potential and may be used for treatment of cardiovascular diseases. Cell biomass grown in bioreactors can be used to produce steroidal saponins. Dioscorea deltoidea plant cell culture contain an extremely complicated mixture of furostanol steroidal saponins. Therefore, it is critical to control the composition of cell cultures for medical use.

Methods

For the extraction of steroidal saponins ultrasound assisted extraction (UAE), multiple successive extraction (MSEM) and reflux extraction (RE) methods were used. LC separation was performed on a reversed phase C18 column in gradient elution mode. MS detection was conducted in the positive electrospray ionization mode (SIM mode).

Results

Optimization of the UAE parameters was performed with use of Latin square experimental design[1]. The following parameters have been optimized: extraction time, organic solvent concentration in extraction solution and the ratio of solvent to sample. Optimal conditions for the extraction of steroidal saponins by the UAE method were: extraction time, 60 min; acetonitrile (water) concentration in extraction solution, 50% (50%); the ratio of solvent to sample, 400 mL/g. As a result, high degree of the extraction of steroidal saponins from the D. deltoidea plant cell suspension culture was achieved. The obtained results were confirmed by MSEM and RE. It was also shown that the UAE method is not suitable for isolation of steroidal glycosides from the D. deltoidea plant material.

Conclusions

For the determination of steroidal saponins from D. deltoidea cell culture comprehensive and representative pre-LC-MS extraction procedure was proposed. Low degree of the steroidal saponins extraction from the plant material is most likely due to the lack of mechanical tissues in the cell culture, in contrast to plant material. The UAE method was tested on a series of 18 cell culture samples of different inoculation date, growing condition and growth period.

Novel Aspect

Comparison of the different approaches performance characteristics have been carried out with respect to finding an optimal pre-LC-MS treatment of *D. deltoidea* cell culture.

References

1. Sarvin B., Fedorova E., Shpigun O., et al., Journal of Chromatography B, 1080, 64-70 (2018).

Acknowledgements

This work was supported by Russian Science Foundation (Grant No. 17-13-01146) for Moscow State University.

For information please contact: stavrianidi.andrey@gmail.com

WP-46 / ANALYSIS OF PESTICIDE RESIDUES IN CANNABIS REGULATED BY CALIFORNIA AND OREGON STATE USING LC/MS/MS WITH DUAL ELECTROSPRAY AND APCI SOURCE

Thomas Becker⁽¹⁾ - *Avinash Dalmia*⁽²⁾ - *Erasmus Cudjoe*⁽³⁾ - *T. Ruthenberg*⁽⁴⁾ - *Josh Ye*⁽³⁾ - *M. Murphy*⁽⁴⁾ - *Feng Qin*⁽³⁾ - *David Welkie*⁽²⁾

PerkinElmer (LAS) Germany GmbH, BLL/FAS-Group, Rodgau-Jügesheim (1) - PerkinElmer Inc., BU, Shelton, CT (2) - PerkinElmer Inc., BU, Woodbridge, ON (3) - SC Labs, Santa Ana, CA (4)

Keywords: Pesticides in Cannabis by LC-MS/MS; Nonpolar and chlorinated pesticides; Dual ESI/APCI source; heated self-cleaning laminar flow interface

Introduction

Over 20 states in the U. S. have legalized the use of recreational or medical cannabis because of therapeutic benefits for ailments such as cancer, multiple sclerosis, and ALS. Since chronic exposure to pesticides in cannabis pose serious health risks, pesticide analysis in cannabis is important. Apart from Oregon state regulatory limits for about 65 pesticides in cannabis, the state of California has issued more stringent regulatory limits for 72 pesticides residues in cannabis flower and edibles. We analyzed all 72 pesticides (including very hydrophobic and chlorinated pesticides generally analyzed by GC/MS) spiked in cannabis flower extracts using LC-MS/MS with ESI and APCI source well below the action limits specified by both Oregon and California State.

Methods

Ground cannabis flower (1 gm) was extracted with acetonitrile (5 mL). The extract was spiked with different levels of 72 pesticides and 10 deuterated internal standards. After that, extract was further diluted 1:1 with acetonitrile containing 0.1 % formic acid. The diluted extract was analyzed directly using LC-MS/MS triple quad with electrospray and APCI source. A 20 min LC-MS/MS method using LX-50 LC, C18 Column and PerkinElmer Qsight 220 MS with electrospray source was used to detect 70 out of 72 pesticides regulated by California state. A 7.5 min LC-MS/MS method using APCI source and C18 column was used to analyze rest of 2 out of 72 pesticides controlled in cannabis by California state.

Results

Normally pesticide analysis in cannabis and other food matrices is done by both GC/MS and LC/MS since some nonpolar and chlorinated pesticides are difficult to ionize with electrospray and APCI ion source used in LC/MS systems. Two different LC-MS/MS methods with electrospray and APCI source were used for low level analysis of all 72 pesticides (including the very hydrophobic and chlorinated pesticides analyzed by GC/MS) in cannabis. The overall sensitivity for most of the pesticides, including those that are normally analyzed by GC/MS was between 1-300 ppb in cannabis, well below regulatory limits set by the state of Oregon and California. The ability to screen and quantitate all 72 pesticides, including the very hydrophobic

and chlorinated GC/MS amenable, in cannabis with LC-MS/MS only with dual ESI/APCI source makes this a novel way of screening and quantitation of pesticides in cannabis and different matrices with a single instrument. The ionization mechanism for analysis of chlorinated pesticides in cannabis using LC-MS/MS system with APCI source would also be presented. The recoveries of majority of the 72 pesticides spiked in grounded cannabis flower were estimated to be between 60-120% with RSD less than 20 %.

Long term stability data for pesticide analysis in cannabis was collected using a triple quadrupole mass spectrometer fitted with dual electrospray ionization source and atmospheric chemical ionization source (APCI) and combined with a heated and self cleaning laminar flow interface. Long term stability data for pesticide analysis in cannabis showed that response RSD over 1 week for majority of pesticides was between 2 to 10 %.

Conclusions

The results demonstrate the successful approach to screen and quantitate pesticides in Cannabis well below the governmental action limits. Very hydrophobic and chlorinated species are also amenable, eliminating the need for further (GC/MS) instrumentation. In addition to this, the maintenance needs of the system are substantially reduced by the heated self-cleaning/laminar flow interface, despite of the presence of challenging matrices such as cannabis.

Novel Aspect

Analysis of GCMS and LCMS amenable pesticides in cannabis using a LC MSMS system with dual ESI/APCI ion source.

WP-47 / IDENTIFICATION OF A NOVEL METABOLITE ORIGINATING FROM JWH-015 USING LCMS-IT-TOF

Özgür Devrim Can⁽¹⁾ - Ümide Demir Özkay⁽¹⁾ - Serkan Levent⁽²⁾ - Nazli Turan Yucel⁽¹⁾ - Fayza Alyu Tekes⁽¹⁾ - Yusuf Ozturk⁽¹⁾

Anadolu University, Faculty of Pharmacy, Department of Pharmacology, Eskisehir, Turkey (1) - Anadolu University, Faculty of Pharmacy, Department of Pharmaceutical Chemistry, Eskisehir, Turkey (2)

Keywords: Narcotic Cannabinoid, JWH-015, Metabolite, LCMS-IT-TOF

Introduction

Despite regulatory efforts; use of “Bonzai” or “Spice” has increased in Turkey due to psychiatric, medical and social features. JWH-015, a narcotic cannabinoid, is one of the active components of Bonzai. Intensive literature survey indicates that studies on the metabolites of JWH-015 are limited to in vitro experiments only [1, 2]. Thus, screening of probable in vivo metabolites of JWH-015 via a non-targeted analysis were aimed in the present study.

Methods

Female Wistar rats of the same age (n=6, ~200-250 g) were used in the experiments. JWH-015 was administered via i.p. route at the dose of 5 mg/kg [2]. Urine samples were collected at varying time points within a 24-hour period and subjected to analysis by using a liquid chromatography mass spectrometry-ion trap-time of flight (LCMS-IT-TOF) instrument (Shimadzu, Japan). Molecular weights of the probable metabolites were calculated and their chemical structures were identified.

Results

In order to determine the pharmacokinetic properties of JWH-015, the analysis of the biotransformation products by LC MS-IT-TOF system were performed. As a result of analyses, a novel in vivo metabolite of JWH-015 was identified for the first time. The novel compound was a glucuronide metabolite of JWH-015, with a chemical name of hydroxy-JWH-015-O-glucuronide. Seven other metabolites with the previously reported chemical structures [1, 2] were also determined in the LCMS-IT-TOF analyses.

Conclusions

Although in vitro metabolism studies give highly accurate results for prediction of possible metabolism pathways and metabolites, in vivo verification of the results is crucial to reveal the real-time reflection of non-living based experiments. Therefore, time dependent determination of novel metabolite, hydroxy-JWH-015-O-glucuronide is very important in terms of pharmacokinetic profile of Bonzai.

Novel Aspect

Discovery of a novel metabolite of JWH-015 gives valuable information about biotransformation pathway of cannabinoid type of compounds.

References

1. Zhang Q., Ma P., Cole R. B., Wang G., Analytical & Bioanalytical Chemistry, 386, 11 (2006).
2. Strano-Rossi S., Anzillotti L., Dragoni S., Pellegrino R. M., Goracci L., Pascali V. L., Cruciani G., Analytical & Bioanalytical Chemistry, 406, 3621-36 (2014).
3. Rizzo V., Carletti F., Gambino G., Schiera G., Cannizzaro C., Ferraro G., Sardo P., Epilepsy Research, 108, 1711-1718 (2014).

WP-48 / SYNTHESIS AND ANTIFUNGAL ACTIVITY EVALUATION OF NEW BENZIMIDAZOLE-TRIAZOLE DERIVATIVES AS POTENTIAL ERGESTEROL BIOSYNTHESIS INHIBITORS USING LC-MS/MS

Nafiz Öncü Can⁽¹⁾ - ***Ulviye Acar Cevik***⁽²⁾ - ***Begüm Nurpelin Sağlık***⁽²⁾ - ***Serkan Levent***⁽²⁾ - ***Yusuf Özkay***⁽²⁾ - ***Zafer Asim Kaplancikli***⁽²⁾

Anadolu University, Faculty of Pharmacy, Department of Analytical Chemistry, Eskisehir, Turkey (1) - Anadolu University, Faculty of Pharmacy, Department of Pharmaceutical Chemistry, Eskisehir, Turkey (2)

Keywords: Ergosterol Quantification, *C. albicans*, Antifungal, LC-MS/MS

Introduction

Ergosterol is the major fungal membrane sterol that regulates membrane fluidity, plasma membrane biogenesis and functions [1]. The inhibition of ergosterol biosynthesis causes a damage in cell membrane of fungi. Hence, quantification of ergosterol level in fungi is essential for development of novel inhibitor candidates [2]. Thus, in the present work, we synthesized new benzimidazole-derivative compounds and developed an LC-MS/MS method for assessment of their effect on ergosterol content of *C. albicans*.

Methods

Sterol extract was prepared from *C. albicans*, incubated with/without inhibitor for 24h [3], and injected (1 µL) to LC-MS/MS system. The mass spectrometric analyses were achieved by employing Nexera XR UFLC system coupled to an LCMS-8040 triple quadrupole mass spectrometer (All from Shimadzu, Japan). MRM was optimized using ergosterol standard, and Labsolutions LC-MS software (Version 5.86, Shimadzu) was used to process the data.

Results

Ergosterol standard was used for quantification of ergosterol in both inhibitor-free (negative control) and inhibitor containing samples. Benzimidazole compounds and reference drugs (ketoconazole and fluconazole) were used at 0.78 µg/mL, 1.56 µg/mL, and, 3.12 µg/mL concentrations. Ergosterol quantity in negative control samples was regarded as 100%. All concentrations were analyzed in quadruplicate, and the results were expressed as mean ± standard deviation. Ergosterol quantification studies indicated that compounds 5i, 5s and reference drugs significantly reduced ergosterol levels at all tested concentrations. Compound 5i displayed 61.74%, 85.41%, and 93.71 % inhibitions at 0.78 µg/mL, 1.56 µg/mL, and, 3.12 µg/mL concentrations, respectively. An inhibitions ratio of 62.05%, 84.48%, and 93.12% were observed for

compound 5s at the same concentrations. Reference drugs indicated similar inhibition potencies to those of compounds 5i and 5s in the range of 63.71% to 94.46%.

Conclusions

We achieved highly acceptable results via developed LC-MS/MS method. Ergosterol quantification studies using *C. albicans* indicated that compounds 5i and 5s and reference agents cause a concentration dependent decrease in the level of ergosterol. Hence, it can be obviously suggested that compounds 5i and 5s have a role in the ergosterol biosynthesis pathway.

Novel Aspect

The developed LC-MS/MS method may hopefully be beneficial for researchers who study on novel antifungal agents as ergosterol biosynthesis inhibitors.

References

1. Zhang Y. Q., Rao R., Virulence, 1, 551-554 (2010).
2. Lupetti A., Danesi R., Campa M., Del Tacca M., Kelly S., Trends in Molecular Medicine, 8, 76-81 (2002).
3. Breivik O. N., Owades J. L., Journal of Agricultural and Food Chemistry, 5, 360–363 (1957).

WP-49 / OPTIMIZED MSⁿ WORKFLOW FOR IMPROVED STRUCTURE ELUCIDATION OF PHARMACEUTICALLY RELEVANT EXTRACTABLES AND LEACHABLES

Seema Sharma⁽¹⁾ - *Kate Comstock*⁽¹⁾ - *Douglas E Kiehl*⁽²⁾ - *Graeme McAlister*⁽¹⁾ - *Ryo Komatsuzaki*⁽¹⁾ - *Caroline Ding*⁽¹⁾ - *Ralf Tautenhahn*⁽¹⁾ - *Derek J Bailey*⁽¹⁾ - *Linda Lin*⁽¹⁾ - *Tim Stratton*⁽¹⁾ - *Shannon Eliuk*⁽¹⁾ - *Iman Mohtashemi*⁽¹⁾ - *Helene Cardasis*⁽¹⁾ - *Jonathan Josephs*⁽¹⁾ - *Vlad Zabrouskov*⁽¹⁾

Thermo Fisher Scientific, San Jose, CA (1) - Eli Lilly and Company, Indianapolis (2)

Introduction

The identification of unknown small molecules, such as impurities, metabolites, degradants, extractables and leachables remains one of the most challenging workflows. Mass spectrometry-based chemical and structural characterization of small molecules greatly benefits from multi stage fragmentation (MSⁿ) coupled with high resolution and high mass accuracy analysis. However, the use of these tools during the untargeted analysis of complex samples is often challenging, owing to the slow precursor interrogation rate and the large number of potential precursors. Here, we describe new data acquisition approaches for the characterization of small molecule extractable and leachable compounds. Building upon sophisticated methods that employ multistage Orbitrap analysis (FTMSⁿ), we dynamically update inclusion and exclusion lists between LC analyses to enable efficient and deep interrogation.

Methods

All experiments were run on a modified Thermo Scientific™ Orbitrap tribrid mass spectrometer. Liquid chromatography separations were carried out on a Thermo Scientific™ Vanquish™ LC system using mobile phases composed of: A: H₂O/0.1% formic acid, and B: ACN/0.1% formic acid with gradient on a Thermo Scientific™ Accucore™ C18 LC column (2.1×100 mm particle size: 2.6µm). A mixture of common plastic and elastomer additives was spiked into matrix and analyzed using LCMS to demonstrate the methodology. This was followed by IPA/water extracts of representative elastomeric components used in biopharmaceutical process. The data was processed using Thermo Scientific™ Mass Frontier™ 8.0 and Thermo Scientific™ Compound Discoverer™ 3.0 software.

Preliminary Data

In this study, a mixture of known additive standards was used to verify the data acquisition features. The IPA/water extract of representative elastomeric components used in pharmaceutical process demonstrates

that using the new features significantly enhances unknown extractable chemical entity identification in a representative solvent extract of an elastomeric component used in pharmaceutical process. The instrument control and data acquisition software for the mass spectrometer was modified to enable automatic generation of inter-run inclusion and exclusion lists based upon the real time component detection in the LC-MS data. Using these inter-run inclusion and exclusion lists, the instrument can selectively trigger MSn scans on features detected in the sample of interest while excluding any distracting background ions. Samples are automatically re-injected to allow for exhaustive MSn data generation for the detected features of interest.

Within each data-dependent LC-MS analysis, we employed sophisticated mass spectrometer methods that are comprised of multiple levels of MSn fragmentation coupled with both high-resolution Orbitrap and low-resolution ion trap analyses. The initial workflow was tested on commercially available standards comprising of plastic and elastomer additives such as antioxidants, slip agents, plasticizers, UV protectants and processing aids. Structure elucidation was conducted based on a newly implemented substructure similarity searching algorithm in Compound Discoverer. Additionally, the MSn data generated for the small molecule compounds of interest was used for generation of a custom spectral library using Mass Frontier 8.0. In the next steps, we will apply this workflow for data analysis of complex unknown small molecule samples.

Novel Aspect

Intelligent data acquisition, processing and interpretation strategy to enable rapid MSn characterization of small molecule unknowns

WP-50 / COMPREHENSIVE METABOLITE IDENTIFICATION OF DM4, A MAYTANSINOID ANTI-TUMOR AGENT, IN S100 FRACTIONS OF HUMAN LIVER MICROSOMES, BY UPLC-QTOF MASS SPECTROMETRY

Elisa Ciuti⁽¹⁾ - *Fabiana Tavani*⁽¹⁾ - *Vera Chirulli*⁽²⁾ - *Chiara Carboni*⁽²⁾ - *Antonio Triolo*⁽¹⁾

Laboratory of Advanced Analytics for Development, Menarini Ricerche Spa, via F. Rismondo 12A, 50131 Firenze (Italy) (1) - Department of Pharmacokinetics and Metabolism, Menarini Ricerche Spa, via . T. Speri 10, 00040 Pomezia (RM, Italy) (2)

Keywords: maytansinoid, anti-tumor, metabolites, microsomes, tandem MS.

Introduction

Maytansinoids are highly cytotoxic molecules; several of them are investigated in cancer therapy, as conjugates with antibodies (ADCs) targeting specific cells. This study describes the identification of the metabolites of the maytansinoid DM4, after incubation with S100 fractions of human liver microsomes (HLM) [1]. Particularly, we wanted to gain better understanding of truly metabolic processes from chemical degradation of the molecule.

Methods

After 1-h incubation, the samples were injected onto a RP C18 UPLC column, and the mass spectrometer was operated in positive ESI alternating low- and high-energy full-scans. Molecular masses of all the DM4-related molecules were obtained by extraction of the masses of typical fragment ions from the high-energy scans (Product Ion Filtration, PIF), and mass defect filtering (MDF) from the low-energy scans. Fragment masses were used for structural elucidation.

Results

Given the reactivity of DM4, the drug-related molecules in HLM incubates were matched with those occurring in control experiments, *i.e.* HLM incubates without the NADPH activator and DM4 incubates in incubation buffer only. The results are as follows.

Incubation in buffer only: DM4 underwent side-chain hydrolysis (*N*-methylglycyl maytansinol), degradation of 1,3-oxazinan-2-one ring of the core structure, SH oxidation to sulphinic and sulphonic acid, hydroxylation of the core ring.

S100 incubation without NADPH: besides the above molecules, the DM4 conjugates with cysteine and, at trace level, with glutathione, were observed.

S100 incubation with NADPH: Sulphinic and sulphonic acid were more abundant. In addition to the above, a metabolite with the side chain SH group replaced by an OH was observed, and a pair of products from Michael reaction of DM4 with 1,4-dihydropyridine-3-carboxamide, resulting from interaction with NADPH or one of its impurities.

Investigations are ongoing to ascertain the nature of the Michael adducts.

Conclusions

The experiments performed on DM4 HLM incubates allowed to distinguish true DM4 metabolites from products arising from merely degradative processes. To the first category belong a sulphinic and sulphonic acid, also partly generated by spontaneous oxidation, cysteine and glutathione conjugation, and side chain SH to OH substitution.

Novel Aspect

Observation of novel DM4 metabolites, first discrimination of P450-dependent processes from other biotransformations and chemical degradation.

References

1. Widdison W. et al., *Molecular Pharmaceutics*, 12, 1762-1773 (2015).

WP-51 / INVESTIGATION OF MONOAMINE OXIDASE INHIBITORY ACTIVITIES OF NEW CHALCONE DERIVATIVES: COMPARISON OF LC-MS/MS AND FLUOROMETRIC METHODS

Ümide Demir Özkay⁽¹⁾ - *Nafiz Öncü Can*⁽²⁾ - *Serkan Levent*⁽³⁾ - *Yusuf Özkay*⁽³⁾ - *Özgür Devrim Can*⁽¹⁾

Anadolu University, Faculty of Pharmacy, Department of Pharmacology, Eskisehir, Turkey (1) - Anadolu University, Faculty of Pharmacy, Department of Analytical Chemistry, Eskisehir, Turkey (2) - Anadolu University, Faculty of Pharmacy, Department of Pharmaceutical Chemistry, Eskisehir, Turkey (3)

Keywords: Monoamine Oxidase, Inhibition, Chalcone, LC-MS/MS

Introduction

Monoamine oxidase (MAO), is a key enzyme that metabolizes monoamine neurotransmitters in peripheral tissues and in the central nervous system [1]. MAO involves two isoforms, *i.e.* MAO-A and MAO-B. MAO-A inhibitors constitute a well-known class of antidepressant drugs, whereas MAO-B inhibitors have importance towards age-related disorders such as Alzheimer's and Parkinson's diseases [2]. In the present work, we applied both fluorometric and LC-MS/MS assays to observe MAO inhibitory potency of new chalcone derivatives.

Methods

MAO enzyme inhibitory effects of newly synthesized chalcone derivatives (2a-2i) [3] were tested in the concentration range of 10^{-3} to 10^{-9} M via a fluorometric method using MAO-A/B Inhibitor Screening Kit (BioVision, CA, USA) according to the manufacturer's guidelines [4]; in addition, a new LC-MS/MS assay method was also developed based on the measurement of tyramine and both methods were compared statistically with each other.

Results

In the fluorometric microplate method, a non-significant inhibition was determined against hMAO-A. In contrast, all compounds indicated higher inhibition profile towards hMAO-B isoform. Compounds 2c, 2e and 2i displayed good inhibition against hMAO-B with lower IC₅₀ values than 1 µM. In order to support the results of enzyme inhibition assay, quantification of substrate (tyramine) level was performed using LC-MS/MS via MRM. Effects of the compounds 2c, 2e and 2i on tyramine level were in agreement with those of compounds observed in fluorometric assay. It was observed that lower inhibitor concentration causes a decrease in the substrate level.

Conclusions

Inhibitors usually react with the enzyme and change it chemically. Hence, in the presence of an inhibitor, the enzyme-substrate complex does not form routinely. Thus, the quantification of tyramine level was performed by an LC-MS/MS method. It was found out that tyramine level increases in the higher concentrations of compounds 2c, 2i and 2f. Consequently, the results of fluorometric kit assay was supported/verified by a more sensitive LC-MS/MS method.

Novel Aspect

Newly reported LC-MS/MS method is very helpful to quantify substrate level in MAO-substrate reaction and to compare results with those of known fluorometric assay.

References

1. Lee J., Lee Y., Park S. J., Lee J., Kim Y. S., Suh Y. G., Lee J., European Journal of Medicinal Chemistry, 130, 365-78 (2017).
2. Youdim M. B. H., Bakhle Y. S., British Journal of Pharmacology, 147, 287-296, (2006).
3. Osmaniye D., Cavuşoğlu K. B., Sağlık B. N., Levent S., Cevik A. U., Atılı O., Ozkay Y., Kaplancıklı Z. A., Molecules, 23, 831 (2018).
4. Fluorometric MAO-A (Catalog No: K796-100) and MAO-B (Catalog No: K796-100) Inhibitor Screening Kit, Data sheet BioVision, USA.

WP-52 / LONG TERM STABILITY OF CANDESARTAN CILEXETIL AFTER MICRONIZATION

Anna Maria Alves de Piloto Fernandes⁽¹⁾ - ***Marcel Menezes***⁽²⁾ - ***Phellipe Honório Amaral***⁽³⁾ - ***Marcos N. Eberlin***⁽¹⁾ - ***Nelci Fenalti Höehr***⁽²⁾

State University of Campinas, Thomson Mass Spectrometry Laboratory, Campinas (1) - University of Campinas, Department of Clinical Pathology, Campinas (2) - University of Campinas, Department of Clinical Pathology, Campinas (3)

Keywords: candesartan, micronization, HRMS, long term stability, degradation products

Introduction

Water solubility is a rate-limiting factor for drug bioavailability. Micronization is the reference technique to increase solubility [1]. In spite of being an old technology, there is little published data about the degradation of pharmaceuticals during this process. Here, ESI-HRMS was used to identify degradation products (DP) of candesartan cilexetil (CC) micronized by milling and stocked for a long time period.

Methods

CC samples was milled for 0.25, 2, 48 or 170 hours. Micronized CC was analyzed 24 months after milling using a Thermo Q-Exactive mass spectrometer at a resolution of 70,000 (@m/z 400). The formation of the DP was evaluated by plotting the ratio of the absolute intensity of the (DP) to that of the CC (IDP/ICC) versus time.

Results

Different times of milling were applied in order to better correlate the micronization process and the formation of the DP of CC. The HRMS analysis of samples 24 months after the micronization process revealed the presence of O-Desethyl Candesartan Cilexetil (Impurity B) and 2-Ethyl-Candesartan Cilexetil (Impurity F). These impurities were identified both in positive (as Na adducts) and negative ion modes. It was verified two different profiles along the long term degradation study: for Impurity B, the ratio IDP/ICC increased seven folds after 48 hours of micronization and remained stable until the end of the experiment. Impurity F, however, showed a continuous growing of the IDP/ICC ratio revealing a close relation between the time of milling during the micronization process and its formation during storage.

Conclusions

Our data indicates a strong correlation between the time of milling and the formation of DP. The kinetic energy transferred to the particles during the process, produces the DP in different extents according to their stability and kinetics of formation. For the first time, the profile of formation of these CC impurities, especially impurity F, were described after a long term storage study.

Novel Aspect

Original HRMS information about the long term stability of candesartan cilexetil after micronization by milling.

References

1. Rasenack N., Muller B. W., Pharm. Dev. Technol., 9, 1 (2004).

WP-53 / MASS SPECTROMETRIC CHARACTERIZATION OF A POLYMERIC DEGRADATION PRODUCT FORMED FROM THE ANTIBIOTIC CEFTAZIDIME

Anna Maria Alves de Piloto Fernandes⁽¹⁾ - *Nilton Soares Camilo*⁽²⁾ - *Elias Paulo Tessaro*⁽¹⁾ - *Marcos N. Eberlin*⁽¹⁾

University of Campinas, Thomson Mass Spectrometry Laboratory, Campinas (1) - University of Campinas, Institute of Chemistry, Campinas (2)

Keywords: ceftazidime, Anti-HIV-1 activity, degradation products, β -lactam antibiotic

Introduction

Degradation products in pharmaceuticals may cause harmful side effects but the new substances generated may have different therapeutic activities. Ceftazidime is a β -lactam antibiotic, which polymerizes during normal storage conditions [1]. Ceftazidime polymer showed a high inhibitory effect on HIV-1 reverse transcriptase [2,3]. In this study mass spectrometry is applied with the objective of revealing the chemical structure of a ceftazidime polymer.

Methods

Ceftazidime (5 g) was heated to 60°C for 210 h in a covered container [2]. A methanolic solution of this polymerized material (1 μ g/mL) was analyzed by direct infusion in a QExactive™ Hybrid Quadrupole-Orbitrap coupled to an HESI II source. MS and MS/MS spectra were acquired at a resolution of 70,000 (@ m/z 400) in the positive ionization mode. MALDI analysis was performed in an Autoflex III Smartbeam MALDI-TOF/TOF instrument using 2,5-DHB as matrix doped with NaCl.

Results

The HRMS showed a distribution of chain lengths typical of a polymeric material, with ions distributed along 300 - 1700 Da mass range. The calculated average molecular weight from the HRMS result was

911.7500 Da. The ions were equally spaced at 74.018 Da units apart. For each oligomeric unit, three different isotopic patterns with mass differences of 5 and 14 Da were observed. MALDI spectra presented the ion clusters between 600 and 1300 Da splitted by a difference of 23 Da suggesting the presence of sodium cationized oligomers.

Conclusions

For the first time, a mass spectrometry based study of ceftazidime polymer has been published. Our mass spectra data, together with other techniques, allows the complete characterization of this interesting degradation product.

Novel Aspect

First mass spectrometry study of a ceftazidime polymeric degradation product.

References

1. Ercanli T., Boyd, D. B., J. Chem. Inf. Model., 46, 1321 (2006)
2. Bcertschi S. W., Cantrel A. S., Kuhfeld M. T., Lorenz U., Boyd D. B., Jcskunos S. R., Antiviral Chem. & Chemother., 8, 353 (1997)
3. Hobi H et al., Antiviral Chem. & Chemother., 12, 109 (2001)

WP-54 / FULL VALIDATION OF METFORMIN IN HUMAN BLOODS USING VOLUMETRIC ABSORPTIVE MICROSAMPLING

Na Ha⁽¹⁾ - ***Kyungho Jang***⁽¹⁾ - ***Song-Hee Han***⁽¹⁾ - ***Seon eui Lee***⁽¹⁾ - ***Min-Gul Kim***⁽²⁾

Center for Clinical Pharmacology and Biomedical Research Institute, Chonbuk National University Hospital, Jeonju (1) - Department of Pharmacology, School of Medicine, Chonbuk National University, Jeonju (2)

Keywords: Microsampling, Volumetric absorptive microsampling, Metformin, LC-MS/MS, Validation

Introduction

Volumetric absorptive microsampling (VAMS) is a novel sampling technique that allows the collection of an accurate volume of blood by dipping a micro-sampler tip. The purpose of this study is to validate the concentration of metformin using LC-MS/MS in VAMS. We demonstrate full validation to ensure that this method is accurate, precise, and stable in metformin.

Methods

A high-performance liquid chromatography with mass spectrometry (LC-MS/MS) was developed and validated for the accurate determination of metformin in human using VAMS. The method was validated for accuracy, precision, sensitivity, recovery, matrix effect, specificity, linearity and stability according to the FDA guideline.

Results

Validation the selectivity of the optimized method was assessed to determine any interfering peak around the retention times of metformin. The calibration curve was calculated by using the linear regression equation of $y=ax+b$, which is the weight ratio of metformin to the peak area of the internal reference material. Calculations of deviations between the measured and theoretical values for the six standard concentrations used in the calibration using the back-calculation method showed that at least 85% of the standard concentrations were within the standard (within 20% of the lowest limit of quantitation, Concentration was within 15%). For standard concentrations that did not meet the acceptance criteria, they were excluded from the calibration curve. The coefficient of determination (r) was higher than metformin 0.9990.

Conclusions

The LC-MS/MS method was used to establish the analytical method for analyzing the concentration of metformin in human blood, and the validity of the established method was verified. In the concentration range of 2 - 100 ng/mL, it showed good linearity and no interference from human blood. The fully validated LC-MS/MS method showed applied to clinical trials such as bioavailability, bioequivalence and pharmacokinetic evaluation of metformin.

Novel Aspect

This is the first paper analyzed metformin by LC-MS/MS using VAMS, a novel technology of microsampling.

References

1. Pinaki S., Uttam B., Animesh G., Amlan K., Bappaditya C., *Chromatographia*, 69, 1243-1250(2009).
2. Srinivasa R., Imran A., Iqbal A., Arindam M., Saral T., *Int J Pharm Pharm Sci*, 5, 463-470 (2013)
3. Philip D., Simon P., Wesley D., Neil S., *J of Pharmaceutical and Biomedical Analysis*, 108, 61-69, (2015)

WP-55 / DETERMINATION OF AN ISONICOTINOHYDRAZIDE DERIVATIVE, A NOVEL POSITIVE INOTROPIC COMPOUND, IN MOUSE PLASMA BY LC-MS/MS AND ITS APPLICATION TO A PHARMACOKINETICS STUDY

Renjie Hui - Jingyi Liu - Bainian Feng

School of Pharmaceutical Science, Jiangnan University, Wuxi

Keywords: isonicotinothiazide derivative, LC-MS/MS, pharmacokinetics, mice

Introduction

Heart failure (HF) is a growing public health problem with high mortality and morbidity [1-2]. The development of HF drugs is a pressing problem [3-5]. In previous works, we presented a new positive inotropic compound, an isonicotinothiazide derivative (AF-HF001), which can directly effect on the myofibril and has a significant effect in the treatment of HF [6-8]. In this study, we employed LC-MS/MS for the qualitative and quantitative analysis of AF-HF001.

Methods

Samples were prepared by one step precipitation with ethyl acetate and stored in acetonitrile. Chromatographic analysis was carried out on a Hypersil Gold C18 column with a gradient mobile phase consisting of acetonitrile and 0.1% aqueous formic acid. The analyte was detected by selective reaction monitoring (SRM) mode with target quantitative ion pair of m/z 292.1→148.2, using praziquantel as the internal standard (IS) m/z 313.1→203.2.

Results

AF-HF001 is a novel positive inotropic compound, which has been reported very recently. The analysis of AF-HF001 in biological matrices is never published. We developed and validated this LC-MS/MS method for the quantitative analysis of AF-HF001 in ICR mice plasma. The one step extraction of ethyl acetate accomplished both drug extraction and protein precipitation, and provided better recovery, stability and simplicity. Then the chromatographic condition with both the composition of mobile phase and the column was also ameliorated to get good sensitivity and selectivity of drug in plasma. With the chosen ion pairs of drug and internal standard, we achieves the quantitative analysis with better precision and sensitivity of AF-HF001 in biological sample. Good linearity ($r=0.995$) was observed over a wide concentration range. The validation of method shows that this bioanalytical method we have developed is sufficient for the analysis of AF-HF001 in mice plasma and will be very essential for the further PK-PD studies of the drug.

Conclusions

A specific and sensitive LC-MS/MS method was developed for the determination of AF-HF001 in mice

plasma for the first time. After a full method validation, it was successfully applied in the preliminary pharmacokinetic studies of AF-HF001 in mice after a single i.p. dose of 45 mg/kg. The analyte was proved as a competitive candidate for the treatment of acute heart failure. Further pre-clinical and clinical pharmacological studies of AF-HF001 are also needed.

Novel Aspect

A full validated LC-MS/MS method of AF-HF001 is described for the first time. It was successfully applied to a very first pharmacokinetics study of AF-HF001 in mice plasma.

References

1. Spirito P, Seidman C, McKenna W, et al. The New England Journal of Medicine 336, 775-785 (1997)
2. Mazurek JA, Jessup M, Heart Failure Clinics 13, 1-19 (2017)
3. Teerlink JR, Metra M, Zaca V, et al. Traditional agents and beyond. Heart Failure Reviews 14, 243-253 (2009)
4. Tamargo J, Amoros I, Barana A, et al. Current Medicinal Chemistry 17, 363-390 (2010)
5. Tamargo J, Lopez-Sendon J, Nature Reviews Drug Discovery 10, 536-555 (2011)
6. Tang C, Xie D, Feng B, Chemical Biology & Drug Design 85, 253-258 (2015)
7. Tang C, Feng B, Journal of China Pharmaceutical University 48, 66-69 (2017)
8. Wang G, Tang C, Yan G, Feng B, Biological and Pharmaceutical Bulletin 39, 207-214 (2016)

WP-56 / LC-MS/MS CHIRAL METHOD FOR SIMULTANEOUS DETERMINATION OF R,S-KETAMINE AND ITS METABOLITES IN RAT PLASMA AND BRAIN SAMPLES FROM PHARMACOKINETIC STUDIES

Monika Janicka - Mikołaj Matłoka - Mateusz Mach - R. Moszczyński-Pętkowski - Ewa Tratkiewicz - Jerzy Pieczykolan - Maciej Wieczorek

Celon Pharma S. A., R&D Department, Łomianki, Poland

Keywords: ketamine, metabolites, chiral chromatography, LC-MS/MS

Introduction

Ketamine is a well-known agent used in anesthesiology since 1960s. The analgesic properties of ketamine were discovered in 1990s. Surprisingly, quite recently this agent also was shown to display unique antidepressive properties. Researchers from Celon Pharma S. A. have conducted comparative analyses of antidepressive effects exerted by R- and S-enantiomers of ketamine and its metabolites in dry powder inhalation as a delivery method.

Methods

Samples (rat plasma and brain) were purified by protein precipitation and centrifugation. Analytes separation was achieved on a chiral Lux® 3 µm Amylose-2 column with a mobile phase consisting of 10 mM ammonium acetate, acetonitrile and 2-propanol. The flow rate was set at 0.5 mL/min. Mass spectrometric detection was performed with use of Agilent 6460 mass spectrometer in positive ionization mode. Ketamine-d₄ was used as an internal standard (IS).

Results

The aim of this study was to develop and validate LC-MS/MS assays for stereoselective quantification of R,S-ketamine and its metabolites R,S-norketamine, (2S,6S)-hydroxynorketamine, (2R,6R)-hydroxynorketamine and R,S-dehydronorketamine in rat plasma and brain samples enabling the characterization of the ketamine pharmacokinetics after administration of a low intravenous (i.v.) or intratracheal (i.t.) dose. Due to very low volume (50 µL) of samples collected during preclinical studies only protein precipitation could be used as a sample preparation. As a result lower limit of quantitation (LLOQ) for all of compounds was determined to be 25 ng/mL. Analytes were separated in one 45 minutes run. Quantitation was

performed using at least two MRM transitions per analyte and IS. The quadratic calibration curves for all analytes were obtained over the concentration range 25–5000 ng/mL. The intra- and inter-day precisions were <7.9% and the accuracy (relative error) was within $\pm 5.3\%$.

Conclusions

The study describes the development and validation of an easy, selective and reliable LC–MS/MS method for simultaneous determination of ketamine enantiomers and its metabolites in low-volume rat plasma and brain samples. The importance of the method is demonstrated by its ability to compare the ketamine metabolites concentrations in these matrices to healthy rats.

Novel Aspect

LC–MS/MS method can differentiate enantiomers of ketamine and its three major metabolites in 45 min run using low volume complex matrix sample.

WP-57 / CHARACTERIZATION OF A NOVEL ZETAMICIN IMPURITY BY LC-MS AND NMR: AN ORTHOGONAL TECHNIQUES APPROACH FOR THE STRUCTURAL ELUCIDATION OF LOW-LEVEL BYPRODUCTS

Elisa Libralesso⁽¹⁾ - *Elisa Ciuti*⁽¹⁾ - *Antonio Triolo*⁽¹⁾ - *Fabio Lamonaca*⁽²⁾ - *Matteo Gentili*⁽³⁾

Menarini Ricerche spa, Advanced Analytics for Development, Firenze (1) - A. Menarini MLS, Analytical Research, Firenze (2) - Giotto Biotech, Organic synthesis, Firenze (3)

Keywords: antibiotics, impurities, LC-MS, NMR

Introduction

Zetamicin is a wide spectrum antibiotic belonging to the aminoglycoside family. Its active principle Netilmicin is made up of a di-aminosaccharide unit (*N*-ethylgaramine) bound to an aminopyran moiety (6-(aminomethyl)-3,4-dihydro-2H-pyran-3-amine). A novel impurity was found in API Zetamicin, while performing a screening for potential related impurities; hence the need to fully characterize it in order to assess its potential toxicity.

Methods

LC-MS analysis provided information on the *m/z* ratio, elemental composition and UV absorption maxima. The isolated impurity was then characterized by mono and bidimensional ¹H NMR, thus refining the structure and assigning the correct isomer among the possible ones. Following the structural determination, the compound was synthetically produced.

Results

Perfect correspondence of the ¹H NMR traces, of MS and MS/MS spectra, and of LC retention times between the process and the synthetic impurities allowed to confirm the assignment to the correct isomer.

Conclusions

This synergic teamwork between orthogonal techniques, such as LC-MS and NMR, was tailored to exploit their unique strengths. In fact, the high sensitivity and separation capabilities of LC-MS allowed to detect and isolate the low-level impurity, and to sketch out of its general structure, while NMR characterization performed on the isolated compound succeeded in the punctual assignment of its structure.

Novel Aspect

Combination of LC-MS and NMR proved effective in the case study, sensibly easing data interpretations, for the strengths of each technique covers the weaknesses of the other.

WP-58 / PREDICTION OF IONISATION EFFICIENCIES IN COMPLICATED BIOLOGICAL MATRICES

Piia Liigand⁽¹⁾ - ***Jaanus Liigand***⁽¹⁾ - ***Filip Cuyckens***⁽²⁾ - ***Rob J. Vreeken***⁽²⁾ - ***Anneli Kruve***⁽³⁾

University of Tartu, Institute of Chemistry, Tartu (1) - Janssen Research and Development, Discovery Sciences, Beerse (2) - Freie Universität Berlin, Institut für Chemie und Biochemie, Berlin (3)

Keywords: ionisation efficiency, electrospray ionisation, metabolites, biological matrices

Introduction

Ionisation efficiency depends on the structure of the compound as well as solvent and setup used. Presently, the relative abundance is evaluated based on the peak area. We demonstrate that the accuracy of such analyses in biological matrices can be significantly improved by predicting ionisation efficiencies.

Methods

All measurements were carried out in ESI negative mode. The logarithm of ionisation efficiencies (logIE) of 10, predominantly pharmaceutical, compounds were predicted in 6 biological matrices (blood, plasma, cerebrospinal fluid, urine, brain and kidney tissue). Charge delocalisation parameter WAPS and degree of ionisation α calculated via COSMO-RS were used as input parameters.

Results

The logIE values in different biological matrices were in good correlation with the logIE values in a solvent. This good correlation suggests that ionisation efficiencies can be predicted in the matrices similarly to the procedure previously proposed for the solvent. To test this, WAPS and α were used to fit the models for predicting logIE values in biological matrices.

The best predictive power of the model was observed in liver tissue homogenate (the root mean square error of the models, sRMSE = 0.49) and the lowest predictive power in urine (sRMSE = 0.80). The correlation between measured and predicted logIE values over all matrices is high, $R^2 = 0.83$ and sRMSE = 0.67 logIE units. Therefore, on average, the mismatch between the predicted and measured ionisation efficiencies is lower than 8 times. This enables semi-quantitative LC/ESI/MS analysis by using a set of calibration compounds to adjust model coefficients for the biological matrix of interest.

Conclusions

Until now equal ionisation efficiencies of all compounds in all matrices are assumed if authentic standards are lacking. For the compounds used in this study, this approach would lead to a mismatch of 660 times with reality. Correcting signals with predicted ionisation efficiencies improves accuracy by almost two orders of magnitude. Additionally, we demonstrate that ionisation efficiencies can be predicted in biological matrices.

Novel Aspect

We present an accurate approach to predict ionisation efficiencies in ESI negative mode for complex biological matrices.

WP-59 / A NOVEL GC-MS METHOD FOR SIMULTANEOUS DETERMINATION OF SEVEN PARABEN DERIVATIVES IN PHARMACEUTICAL AND COSMETIC PRODUCTS CONSUMED IN TURKEY

Saniye Ozcan⁽¹⁾ - ***Serkan Levent***⁽²⁾ - ***Aysun Geven***⁽¹⁾ - ***Nafiz Öncü Can***⁽¹⁾

Anadolu University, Faculty of Pharmacy, Department of Analytical Chemistry, Doping and Narcotic Compounds Analysis Laboratory, Eskisehir (1) - Anadolu University, Faculty of Pharmacy, Department of Pharmaceutical Chemistry, Doping and Narcotic Compounds Analysis Laboratory, Eskisehir (2)

Keywords: GC-MS, parabens, isomers, pharmaceuticals, cosmetics

Introduction

Due to official regulations on the use of Alkyl esters of p-hydroxybenzoic acid (parabens) [1], the presented study was developed for quantitative determination of methylparaben (MP), ethylparaben (EP), n-propylparaben (NPP), i-propylparaben (IPP), n-butylparaben (NBP), i-butyl paraben (IBP) and benzylparaben (BP) in some pharmaceutical and cosmetic products consumed in Turkey using gas chromatography-mass spectrometry (GC-MS) - isotope dilution method.

Methods

GC was carried out using a highly inert 95% dimethylpolysiloxane, 5% diethylpolysiloxane phase column (Rxi®-5ms, 30 m × 0.25 mm × 0.25 µm film thickness, Restek, USA). The ion trap mass spectrometer was operated in the electron impact - positive ionization mode (+70 eV) using an external ionization configuration. In the full scan mode the mass range was varied from 35 to 500 *m/z* at 0.6 s scan⁻¹.

Results

Analyte concentrations were determined by isotope dilution method using ¹³C¹³-MP-EP-PP-BP solutions as an internal standard to quantitative the unlabeled paraben derivatives. The obtained results for linearity, precision, and accuracy of the proposed method support their potential use for simultaneous identification and quantitative determination of the MP, EP, NPP, IPP, NBP, IBP and BP in samples. The method was validated according to ICH Q2(R)1 regulations [2], and system suitability parameters were checked according to the USP.

The amount of each of paraben was calculated using the calibration equations in 17 pharmaceutical and cosmetic products. Among the samples, none of the above-mentioned parabens were found in 11 of the analyzed products; both MP and NPP were detected in 4 of the products.

Conclusions

A simple, selective and highly accurate isotope dilution - GC-MS method was developed and applied to analyze the MP, EP, NPP, IPP, NBP, IBP and BP in cosmetics and pharmaceutical products. Further studies with more products of different compositions and different manufacturers have to be conducted to make a better assessment overall possible hazard from originating from pharmaceutical and cosmetic products.

Novel Aspect

To the best of our knowledge, this is the first study in which pharmaceutical or cosmetics products from Turkish market were analyzed for parabens content using GC-MS.

References

1. Moos R. K., Koch H. M., Angerer J., Apel P., Schroter-Kermani C., Bruning T., Gehring M. K., International Journal of Hygiene and Environmental Health, 218, 666-74 (2012).
2. ICH Harmonised Tripartite Guideline: Validation of Analytical Procedure Text and Methodology Q2 (R1), International Council on Harmonisation of Technical Requirements for Registration of Pharmaceuticals for Human Use, (2005).

WP-60 / ANTICANCER ACTIVITY EVALUATION OF NEW PHORTRESS ANALOGUES AND IDENTIFICATION OF BIOACTIVE METABOLISM PRODUCTS USING LCMS-IT-TOF

Yusuf Özkay⁽¹⁾ - **Derya Osmaniye**⁽¹⁾ - **Begüm Nurpelin Sağlık**⁽¹⁾ - **Serkan Levent**⁽¹⁾ - **Sinem Ilgin**⁽²⁾ - **Zafer Asim Kaplancikli**⁽¹⁾

Anadolu University, Faculty of Pharmacy, Department of Pharmaceutical Chemistry, Eskisehir, Turkey (1) - Anadolu University, Faculty of Pharmacy, Department of Pharmaceutical Toxicology, Eskisehir, Turkey (2)

Keywords: Phortress, Benzoxazole, Anticancer, Active Metabolite, LCMS-IT-TOF

Introduction

Phortress, including benzothiazole moiety, is a prodrug displaying anticancer activity. The active metabolite of phortress is a potent agonist of the aryl hydrocarbon receptor and switches on cytochrome P450 CYP1A1 gene expression [1]. In the present study, new analogues of phortress were synthesized and investigated for anticancer activity. In order to identify a probable active metabolite, biotransformation studies were conducted and results were examined by mass spectroscopic analyses.

Methods

Target benzoxazole derivatives (3a-3p) were synthesized [2], and anticancer activity tests were performed on colon (HT-29), breast (MCF7), lung (A549), liver (HepG2) and brain (C6) carcinoma cell types [3]. Induction potential of the compounds 3m and 3n on CYP1A1/2 enzymes was additionally investigated. Biotransformation studies for 3m and 3n were evaluated by LCMS-IT-TOF instrument (Shimadzu, Japan). Moreover, molecular docking studies for compound 3n and its some metabolites were performed.

Results

Anticancer potentials of phortress analogs 3a-3p were assessed by cytotoxicity, DNA synthesis inhibition, and apoptosis induction assays. Compounds 3m and 3n were the most cytotoxic agents, which inhibited DNA synthesis of carcinogenic cell lines, significantly. It was also determined that these compounds have an effect to induce apoptotic cell death. In the biotransformation studies compounds 3m and 3n were incubated with HepG2 cell line for 24 h. Metabolite identification studies were performed using LCMS-IT-TOF system. Proposed biotransformation pathway displayed that there are metabolites (3mM1, 3mM2, 3nM1, and 3nM1), which undergo the same metabolic reactions as active metabolites of phortress. CYP1A1/2 enzyme induction and molecular docking studies gave similar results to those of phortress.

Conclusions

As synthesized compounds thought to be prodrugs, their metabolic pathways were evaluated by applying biotransformation studies. Identification of very similar metabolites to active metabolites of phortress recommended that 3m and 3n have same mechanism of action with that of phortress. Furthermore, results of CYP1A1/2 enzyme induction and molecular docking studies also supported this suggestion.

Novel Aspect

Two novel anticancer drug candidates, which have stronger activity than phortress, were developed along with their active metabolites determined by LCMS-IT-TOF analyses.

References

1. Bradshaw T. D., Stone E. L., Trapani V., Leong C. O., Matthews C. S., tePoele R., Stevens M. F. G., Breast Cancer Research and Treatment, 110, 57-68 (2008).
2. Aiello S., Wells G., Stone E. R., Kadri H., Bazzi R., Bell D. R., Stevens M. F. G., Matthews C. S., Bradshaw T. D., Westwell A. D., Journal of Medicinal Chemistry, 51, 5135-5139. (2008).
3. Osmaniye D., Levent S., Karaduman A. B., Ilgin S., Ozkay Y., Kaplancikli Z. A., Molecules, 23, 1054 (2018).

WP-61 / INVESTIGATION OF ISOMERIC METABOLITE STRUCTURES USING A COMBINATION OF IMS, MS/MS AND CHEMOMETRICS: APPLICATION TO BIOTRANSFORMATION.

Cécile Palaric⁽¹⁾ - *Sébastien Rigaud*⁽¹⁾ - *David Mathiron*⁽¹⁾ - *Roland Molinié*⁽²⁾ - *Jean-Yves Beaumal*⁽³⁾ - *Tristan Renaud*⁽⁴⁾ - *Serge Pilard*⁽¹⁾

University of Picardie Jules Verne (UPJV), Analytical Platform (PFA), Amiens (1) - University of Picardie Jules Verne (UPJV), BIOPI-EA3900, Amiens (2) - Technologie Servier, Synthesis, Orléans (3) - Technologie Servier, A. S. M. P., Orléans (4)

Keywords: Biotransformation; Position of substitution; MS/MS; Ion mobility; Chemometrics.

Introduction

Biotransformation is used to build molecules, similar to active principal ingredient, in the drug discovery program [1]. The structural identification of metabolites remains challenging when isomeric substituted molecules are involved. In this work, we present a methodology which combines, travelling-wave ion mobility separations (TWIMS) [2], MS/MS fingerprints and chemometrics [3], to successfully assign a number of isomeric hydroxylated metabolites.

Methods

A model molecule bearing an aza-bicyclo moiety and 6 of its hydroxylated metabolites, previously isolated, were used for our methodology development. Retention times (tR) from UPLC/HRMS separation, instrument measured drift times (tD) from TWIM separation of the protonated molecules, accurate m/z and intensity measurements of fragment ions from MS/MS spectra, were collected.

Results

The tR, drift times and MS/MS data set of our standard molecules led to coherent information regarding hydroxylated metabolites substitution position, shape, polarity and stereochemistry. Statistical processing of MS/MS fragment ions intensities, using multivariate analysis (SIMCA-P), enabled a fast and unambiguous assignment of the metabolites structures, which could allow the transposition of our approach to investigate complex bioconversion mixtures.

Conclusions

To validate our workflow, the UPLC/HRMS screening of 12 fungi stains followed by a principal component analysis (PCA) were performed. After the selection of the best candidate, the hydroxylated metabolites content was quickly elucidated with the previously described tools.

Novel Aspect

Investigation of isomeric compounds using a combination of ion mobility, MS/MS and chemometrics.

References

1. Kebamo S., Tesema S., Geleta B., Journal of Drug Metabolism & Toxicology, 6:5 (2015).
2. Dear G. J., Munoz-Muriedas J., Beaumont C., Roberts A., Kirk J., Williams J. P., Campuzano I., Rapid Communications in Mass Spectrometry, 24, 3157–3162 (2010).
3. Otto M., Chemometrics, Wiley-VCH, Third Edition (2017).

WP-62 / COMPLEX APPROACH FOR ISOLATION AND HPLC-MS DETERMINATION OF STEROIDAL GLYCOSIDES FROM PLANT MATERIAL AND SUSPENSION CELL CULTURE

Igor Rodin - Boris Sarvin - Andrey Stavrianidi - Timur Baygildiev

Lomonosov Moscow State University, Department of Chemistry, Moscow

Keywords: extraction, HPLC-MS, steroidal saponins, experimental design

Introduction

For thousands of years, herbal medicines made from terrestrial plants were used by humans to prevent and treat diseases. Extraction is the crucial first step in the analysis of medical plants. Researches have shown that most of the therapeutic properties of *D. deltoidea* and *T. terrestris* are due to the presence of steroidal glycosides. The main objectives of this study were to develop a complete strategy of isolation and HPLC-MS determination of steroidal glycosides.

Methods

The HPLC separation was conducted on a C18 column. The separation was carried out in a gradient elution mode. Detection were performed by MS in positive ion (SIM mode). Ultrasound-assisted (UAE), refluxing (RE), low pressure refluxing (LPRE) and Soxhlet (SE) extraction techniques were developed for determination steroidal saponin content. The results were confirmed using the multiple successive extraction method and spiking with standards before extraction.

Results

Extraction parameters for UAE, RE, LPRE and SE methods were optimized by single-factor-experiment for isolation of steroidal glycosides from *T. terrestris* [1]. Extraction parameters for UAE method were optimized by Latin Square experimental design at 4 levels for 4 factors for isolation of steroidal glycosides from *D. deltoidea* plant material and cell culture [2]. Thermal decomposition of protodioscin was observed after long term high temperature extraction process. The accuracy for all studied extraction techniques was confirmed by spiking of the plant material with standards before extraction and HPLC-MS analysis and the analytical method was validated for linearity, limits of detection, limit of quantification, precision and accuracy. Also, a high degree of extraction was confirmed by multiple successive extraction method (MSEM) and RE methods. At the same time, it was shown that the UAE method is completely inappropriate for a similar task when working with plant material. This may be due to the lack of mechanical tissues in the cell culture, in contrast to plant material.

Conclusions

In the course of the work cycle, a study was made of the patterns and features of extraction, chromatographic separation followed by mass spectrometric determination of the steroidal glycosides from plant material and cell cultures. The impossibility of using the UAE extraction method for working with plant material was also shown.

Novel Aspect

The method for HPLC-MS determination of steroid glycosides was developed and validated. The optimization and comparison of extraction techniques were performed.

References

1. Sarvin B., Stekolshchikova E., Rodin I., Stavrianidi A., Shpigun O., JARMAP, 8, 75-82 (2018).
2. Sarvin B., Fedorova E., Shpigun O., Titova M., Nikitin M., Kochkin D., Rodin I., Stavrianidi A., Journal of Chromatography B, 1080, 64-70 (2018).

Acknowledgements

This work was supported by Russian Science Foundation (Grant No. 17-13-01146) for Moscow State University.

For information please contact: igorrodin@yandex.ru

WP-63 / USE OF HIGH RESOLUTION MASS SPECTROMETRY FOR TROPANE ALKALOIDS DETERMINATION AND UNTARGETED ANALYSIS IN DRUGS AND HOMEOPATHIC PRODUCTS FROM ATROPA BELLADONNA SEED EXTRACTS

Roberto Romero González - Jesús Marín Sáez - Antonia Garrido Frenich

Research Group "Analytical Chemistry of Contaminants", Department of Chemistry and Physics, Research Centre for Agricultural and Food Biotechnology (BITAL), University of Almeria, Agrifood Campus of International Excellence, ceiA3, E-04120 Almeria, Spain

Keywords: Homeopathic products, tropane alkaloids, *Atropa belladonna*, Orbitrap-MS, untargeted analysis.

Introduction

Tropane alkaloids are compounds with therapeutic properties naturally obtained from Solanaceae plants as *Atropa Belladonna*, which has been used in drugs and homeopathic products. The main concern related to homeopathic products is the scarce information about their composition as well as recent cases of secondary effects [1]. New methods are necessary to control these products, including untargeted analysis.

Methods

The extraction method was based on solid phase extraction (SPE) for homeopathic products and then, the recomposed extract was analyzed by liquid chromatography coupled to high resolution mass spectrometry (Exacte-Orbitrap). For drugs (ointment) and *Atropa Belladonna* seeds only dilution was used prior injection. A homemade database with 41 tropane alkaloids was created and applied. Compound Discover 2.1 was used to search untargeted compounds.

Results

The method was optimized using a SPE procedure to concentrate the compounds present in homeopathic products, as well as the dilution step for drugs and *Atropa Belladonna* seeds. The method was validated evaluating parameters as specificity, linearity, accuracy, precision and limits of detection (LOD) and quantification (LOQ). The proposed method allows the determination of a large range of tropane alkaloids in *Belladonna* seeds and ointment from *Belladonna* extract at concentration up to 337.5 mg/kg for atropine, while in homeopathic products only atropine, scopolamine, tropine and apoatropine were detected. When untargeted analysis was performed, other tropane alkaloids (benzoylecgonine) found in related plants [2] were also detected. In addition to tropane alkaloids, aminoacids, considered as precursors of tropane alkaloids (phenylalanine and ornithine), were confirmed in the seeds as well as sugars (glucose, fructose and lactose) in the homeopathic drugs.

Conclusions

The developed method provides the quantification of a large range of tropane alkaloids. High resolution mass spectrometry allows the retrospective analysis of the samples and the available software makes possible untargeted analysis, detecting other tropane alkaloids, precursors as aminoacids and other compounds present in the drugs.

Novel Aspect

A suitable methodology for the comprehensive characterization of homeopathic products from *Atropa Belladonna* applying targeted and nontargeted approaches is provided.

References

1. Glatstein, M., Danino, D., Wolyniez, I., Scolnik, D., *American Journal of Therapeutics*, 21, e196–e198 (2014).
2. Jirschitzka, J., Schmidt, G. W., Reichelt, M., Schneider, B., Gershenzon, J., D'Auria, J. C., *Proceedings of the National Academy of Sciences*, 109, 10304–10309 (2012).

WP-64 / ALL-IN-ONE METHOD FOR PROFILING AND QUANTIFICATION OF POLYSORBATE, BY COMBINED DATA-INDEPENDENT AND DATA-DEPENDENT TANDEM MS

Serena Staccioli - Elisa Ciuti - Antonio Triolo

Laboratory of Advanced Analytics for Development, Menarini Ricerche, via Francesco Rismondo 12A, 50131 Firenze (Italy)

Keywords: Polysorbate, quantification, characterization, tandem MS, all-in-one method

Introduction

Polysorbate (PS) is a polyoxyethylene (POE) type non-ionic surfactant derived from esterification of sorbitan, made from a very complicated mixture of polymers, widely used as excipient.

Analytical challenges are: structural heterogeneity caused by different sorbitan core structures esterified in different positions with a mix of fatty acids; small amount in formulations, potential interferences from excess of drug substance and excipients; no UV absorbance.

Methods

All the information is obtained in one HPLC-MS run using a Q-ToF instrument. The method includes a low-energy survey MS scan, to obtain information on the intact components and the molecular weight distributions of the involved species; a data-independent MS/MS scan to detect characteristic fatty acids marker ions; a data-dependent MS/MS of selected precursor ions, triggered by detection of a given fatty acid marker in the data-independent scan.

Results

We developed and validated an analytical method for simultaneous profiling and quantification of polysorbate 80 (having oleic acid as the main fatty acid component) either as neat excipient or contained in pharmaceutical and biopharmaceutical formulations.

The method provides identification of components, their apparent percentages, the fatty acids profile, the absolute quantification of this excipient in formulations, the detection and identification of its degradants.

The developed method finds different application in quality control, comparison of suppliers, monitoring batch-to-batch consistency, regulatory documentation, understanding preformulation and stability issues.

Conclusions

The developed method is suitable for rapid characterization and quantification of neat polysorbate and polysorbate in pharma and biopharma formulations, yielding simultaneously components identity, their apparent percentages, the fatty acids profile and its absolute quantification in different kinds of solid and liquid formulations.

Novel Aspect

First description of a suitable qualitative/quantitative method for comprehensive characterization of PS.

For information please contact: sstaccioli@menarini-ricerche.it

WP-65 / RAPID HPLC-MS/MS EVALUATION OF THE BATCH-TO-BATCH CONSISTENCY OF THE PAYLOADS SUBSTITUTION OF A LYSINE ANTIBODY-DRUG-CONJUGATES (ADCS).

Fabiana Tavani - Antonio Triolo

Menarini Ricerche SpA, Advanced Analytics for Development, Florence

Keywords: ADC; Conjugation site analysis; Lysine conjugation; Data-independent acquisition; Batch-to-batch consistency

Introduction

The aim of this work is the determination of the drug substitution sites of a lysine antibody-drug conjugate (ADC) by a very specific LC-MS analysis of its peptide digests [1][2]. Conjugated peptides were detected with high sensitivity and specificity by monitoring the characteristic reporter fragment ions of the cytotoxic molecule. Moreover, the batch-to-batch consistency of the conjugation was studied to evaluate the reproducibility of the production process.

Methods

The ADC was denatured and digested with a Trypsin/Lys-C mix (E/P ratio 1:25). The instrument was a Synapt G2 HDMS mass spectrometer interfaced with a Waters Acquity UPLC H-Class System, using an Acquity Peptide CSH UPLC column. The data acquisition was in MSE mode. Conjugated peptides were detected by extracting the payload fragment ions from the high energy trace. The results were elaborated by UNIFI software (Waters).

Results

From the automatic data processing with UNIFI software, about 30 payload-substituted peptide were detected, which were then visually inspected to confirm their identity. The conjugated peptides were characterized by their accurate masses, retention times and by their MS/MS sequence, allowing unambiguous detection of the substituted residues [3].

Subsequently, eight batches of ADC were treated with the same experimental method and the extracted ion chromatograms of the payload reporter ions were compared: in this way it was possible to highlight similarities and differences among the obtained chromatographic profiles, and hence the similarity degree of the substitution sites among all batches. The evaluation of the obtained chromatograms showed that they were practically superimposable, demonstrating excellent batch-to-batch consistency of the payload substitution.

Conclusions

The payload-substituted peptides in the digestion mixture were detected with high specificity, allowing substitution site determination based on sequence information from the MS/MS data. Moreover, the comparison of the profiles of the extracted chromatograms of the payload fragment ions allows to evaluate the reproducibility of the production process and to identify any differences, in the substitution sites, among several batches.

Novel Aspect

Description of a quick and simple method to assess batch-to-batch reproducibility of the conjugation process in an ADC.

References

1. L. Chen et al, MABS, 8: 1021 (2016).
2. Q. Luo et al, Analytical Chemistry, 88: 695 (2016).
3. H. W. Chih et al, J. Pharm. Sci, 200: 2518 (2011).

WP-66 / MALDI TOF MASS SPECTROMETRY FOR THE CHARACTERIZATION OF ANTITUMOUR RUTHENIUM(II) POLYPYRIDYL COMPLEXES

Romana Masnikosa⁽¹⁾ - **Filip Veljkovic**⁽¹⁾ - **Boris Rajcic**⁽¹⁾ - **Ana Rilak**⁽²⁾ - **Suzana Velickovic**⁽¹⁾

University of Belgrade, Vinča Institute of Nuclear Sciences, PO Box 522, 11001, Belgrade, Serbia (1) - University of Kragujevac, Faculty of Science, R. Domanovica 12, P. O. Box 60, 34000 Kragujevac, Serbia (2)

Keywords: MALDI-TOF-MS, ruthenium(II) complexes, instrumental parameters

Introduction

Ruthenium(II) polypyridyl compounds, such as Ru(II) terpyridine complexes, are among lead compounds in anticancer drug discovery [1]. Their characterization after synthesis is, thus, necessary, but also presents a challenge. To this end we used matrix-assisted laser desorption/ionization time of flight mass spectrometry (MALDI-TOF-MS). Our main goal was to find how instrumental parameters affect the mass spectra of $[\text{Ru}(\text{Cl-Ph-terpy})(\text{N-N})\text{Cl}]^+$ complex ions.

Methods

Experiments were performed on the MALDI TOF Voyager-DE PRO (Sciex, USA) instrument, equipped with a delayed ion-extraction device and a pulsed nitrogen laser (337 nm 3 ns pulse width) operating at 20.0 Hz laser repetition rate. Solution of dithranol in tetrahydrofuran was employed as a matrix. A 0.5 μL of the matrix was placed onto a stainless steel target plate and laid over by 0.5 μL of a sample solution.

Results

MALDI MS has been successfully applied for identification and characterization of transition metal compounds and organometallics [2-5]. In spite of this, data regarding the effect of matrix, laser intensity and analysis mode on mass spectra of transition metal complexes is limited [6]. In the study presented herein we optimised the values of input variables: the laser intensity, the extraction delay time, the accelerating voltage, the grid voltage and the number of laser shots for the characterization of the three complexes. The complexes have a general formula $\text{mer}[\text{Ru}(4'\text{-chlorophenyl-}2,2':6',2''\text{-terpyridine})(\text{N-N})\text{Cl}]\text{Cl}$, where N-N is a bidentate ligand: 1) 1,2-diaminoethane (en), 2) 1,2-diaminocyclohexane (dach) or 3) 2,2'-bipyridine (bpy). The fact that Ru possesses seven stable isotopes aggravated the analysis of mass spectra of its compounds. Since the studied complexes contain an inner chloride, the isotopic distribution of Cl must also be taken into an account.

Conclusions

In the mass spectra of the studied Ru(II) complexes the peaks of molecular (non fragmented) ions were visible only upon the careful selection of the instrumental parameters. Otherwise fragment ions dominated the mass spectra, which precluded the determination of the precise molecular masses of the complexes. The accelerating voltage, the grid voltage and the extraction delay time had the greatest influence on the mass spectra and isotopic resolution of the peaks.

Novel Aspect

Tuning of instrumental parameters of the MALDI TOF mass spectrometer for characterization of novel anticancer Ru(II) terpyridine compounds.

References

1. Lazarevic T., Rilak A., Bugarcic D., Eur. J. Med. Chem., 142, 8(2017)
2. Dopke C., Treichel M., Vestling M., Inorg. Chem., 37, 1272(1998).
3. Hunsucker W., Watson C., Tissue M., Rapid Commun. Mass Spectrom., 15, 1334 (2001).
4. Grishin I, J. Anal. Chem., 68, 1123 (2013).
5. Baranov E., Fukin G., Balashova T., Pushkarev A., Grishinb I., Bochkarev M., Dalton Trans., 42, 15699 (2013).
6. Wyatt M., J. Mass Spec., 46, 712 (2011).

WP-67 / ENERGY RESOLVED TANDEM MASS SPECTROMETRY EXPERIMENTS FOR RESOLUTION OF ISOBARIC COMPOUNDS: A CASE OF CIS/TRANS ISOMERISM

Laura Braconi - Marta Menicatti - Silvia Dei - Elisabetta Teodori - Gianluca Bartolucci

Università degli Studi di Firenze, Department of Neurofarba, Sesto Fiorentino, Italy

Keywords: linear equations of deconvolution analysis, LEDA, energy resolved mass spectrometry, isomers characterization, multi-drug resistance

Introduction

A series of cis/trans ester isomers, that showed high efficacy to reverse the acquired resistance of cancer cells during chemotherapeutic therapy was studied as candidate drugs [1]. In order to evaluate their stability in human plasma, a liquid chromatography tandem mass spectrometry (LC-MS/MS) method was

developed together with a mathematical algorithm, that is able to identify the isomers present in analyzed samples, without their chromatographic separation.

Methods

The LC-MS/MS analyses were carried out using a triple quadrupole system equipped with an electrospray source (ESI), operating in positive ion mode. The employed chromatographic parameters were finely tuned to minimize the run time. Therefore, a short C18 Column (3 cm) and fast elution gradient were used. Furthermore, a matrix of linear equations (LEDA), based on abundance ratios of product vs precursor ions, was applied to allow the isomers distinction[2-3].

Results

Using the reported conditions, a co-elution of the studied isomers was observed. Analysis of isobaric molecules, and especially isomers, by tandem mass spectrometry is often complicated by the similarity between their fragmentation patterns and it is common that the same MS/MS product ions are present in the spectra of all the isomers. Hence, to find some differences among cis/trans isomers, a series of experiments based on energy resolved MS/MS (ERMS) was carried out. Knowing the characteristic ratios of pure isomer, LEDA algorithm elaborates each MS/MS acquired spectra, assigning the correct abundance to the identified isomer. Therefore, in case of an isomers' mixture, LEDA is able to resolve the co-eluting peak determining the relative concentration of each isomer eventually present.

The LEDA reliability was checked during stability studies. The results proved the effectiveness of correct assignment of the isomer present in the sample. This ability was essential to verify the behavior of the studied isomers during stability tests.

Conclusions

The proposed approach, combining a LC-MS/MS method with a post-processing mathematical algorithm (LEDA), was able to deconvolute a mixture of isomers from chromatographically unresolved peaks. The results demonstrated that this combination was suitable for the evaluation of isomers' degradation profiles in human plasma samples. Indeed, the stability experiments have highlighted a different behavior between cis/trans isomers.

Novel Aspect

The LEDA approach has the advantage that isomers can be quantified without the need of LC separation or additional specialized ion mobility instrumentation.

References

1. Teodori E., Dei S., Coronello M., Floriddia E., Bartolucci G., Manetti D., Romanelli M. N., Porqueras D. S. D., Salerno M., Eur. J. Med. Chem. 127, 586-598 (2017).
2. Menicatti M., Guandalini L., Dei S., Floriddia E., Teodori E., Traldi P., Bartolucci G., Rapid Commun. Mass Spectrom. 30, 423-432 (2016).
3. Menicatti M., Guandalini L., Dei S., Floriddia E., Teodori E., Traldi P., Bartolucci G., Eur. J. Mass Spectrom. 22, 235-243 (2016).

WP-68 / STABILITY STUDY IN PLASMA OF NOVEL ESTER-LINKED NSAID - CARBONIC ANHYDRASE INHIBITOR HYBRIDS

Silvia Bua - Marta Menicatti - Paola Gratteri - Alessio Nocentini - Claudiu T. Supuran - Gianluca Bartolucci

University of Florence, NEUROFARBA Dept., Pharmaceutical and Nutraceutical section, via Ugo Schiff 6, 50019 Sesto Fiorentino (Florence), Italy

Keywords: Stability, NSAID, Carbonic Anhydrase, Hybrids

Introduction

We reported a series of amide-based hybrid compounds incorporating 6- and 7-substituted coumarins and some NSAIDs as agents for the management of RA. [1] These derivatives were proved to exert remarkable anti-inflammatory action in spite of non-significant amide-bond cleavage in vivo. [1-3]

As a rational development, we report herein all such derivatives ester bioisosters to investigate the contribution of a more cleavable linker in terms of plasma stability.

Methods

The kinetic studies were carried out through a stopped-flow assay. [4] The LC-MS/MS analysis were carried out using a triple quadrupole system equipped with an electrospray source (ESI), operating in positive-ion mode. The employed chromatographic parameters were finely tuned to minimize the run time. Furthermore, a matrix of linear equations (LEDA), based on abundance ratios of product vs precursor ions, was applied to allow the isomers distinction [3].

Results

The ester hybrids were shown to inhibit the target human CA isoforms more efficiently than the corresponding amides. Using the LC reported conditions, a co-elution of the studied isomers were observed. Knowing the characteristic ratios of pure isomer, LEDA algorithm elaborates each MS/MS spectra acquired, assigning the correct abundance to the identified isomer. Therefore, in case of a isomers mixture, LEDA is able to resolve the co-eluting peak determining the relative concentration of each isomer eventually present. The LEDA reliability was checked during stability studies. The results proved the effectiveness of correct assignment of the isomer present in the sample. The drug stability studies demonstrate that the amides derivatives resulted to be stable in both human and rat plasma. Likewise, most esters derivatives exhibited human plasma stability, with the exception of Ketoprofen and Diclofenac derivatives, showing that the ester bond can undergo hydrolytic cleavage in human plasma depending on the NSAID fragment.

Conclusions

Docking studies displayed the capability of the overall hybrid to bind to COX-1 and COX-2. Unlike amides, the ester bond can undergo hydrolytic cleavage in human plasma depending on the NSAID fragment. In vitro functional investigations and in vivo studies are currently ongoing to finally determine the mechanism of action of such derivatives.

The proposed LC-MS/MS approach with LEDA tool, was able to separate the isomers eventually present from LC unresolved peaks.

Novel Aspect

Introduction of ester/amide-linked NSAIDs-coumarins hybrids for the management of RA. The use of LEDA tool to isomers analysis without LC or ion mobility separation.

References

1. Bua S., Di Cesare Mannelli L., Vullo D., Ghelardini C., Bartolucci G., Scozzafava A., Supuran C. T., Carta F., *J Med Chem*, 60, 1159-1170, (2017)
2. Neri D., Supuran C. T., *Nature Rev Drug Disc*, 10, 767-777, (2011)
3. Menicatti M., Pallecchi M., Bua S., Vullo D., Di Cesare Mannelli L., Ghelardini C., Carta F., Supuran C. T., Bartolucci G., *J Enzyme Inhib Med Chem*, 33, 671-679, (2018)
4. Khalifah RG, *J Biol Chem*, 246, 2561-73, (1971)

WP-69 / IDENTIFICATION OF DIPROPYLAMINOPRETADALAFIL USING HIGH RESOLUTION MASS SPECTROMETRY UNDER MULTI-STAGE FRAGMENTATION PATHWAYS

Agnes Chin - Chee Leong Kee

Keywords: Dipropylaminopretadalafil, Diethylaminopretadalafil, HRMS, multistage fragmentation

Introduction

Phosphodiesterase type 5 (PDE-5) inhibitors and their analogues are commonly detected in illegal health products used to treat erectile dysfunction in men. The emerging health risks have become prominent with the first fatal case associated with intoxication by desmethylcarbodenafil reported in 2017 [1]. Therefore, it is crucial that the laboratory is equipped with the necessary equipment and knowledge to structural elucidate the structure of new analogues. In this study, a new tadalafil analogue was detected in a routine analysis in the laboratory. The structure of this unknown was compared to that of a known tadalafil analogue, diethylaminopretadalafil, using liquid chromatography mass spectroscopy (LC-MS) and high-resolution mass spectroscopy (HRMS).

Methods

The unknown PDE-5 inhibitor analogue was isolated by liquid-liquid extraction and purified by recrystallization process. The isolated compound was then dissolved in acetonitrile for Liquid Chromatography-Mass spectroscopy (LC-MS). The sample solution was further diluted with 0.1% formic acid in water: 0.1% formic acid in acetonitrile = 4:1 for direct infusion using a Thermo Fischer Scientific LTQ Orbitrap XL FTMS. The infusion flow rate was 10 μ L per minute. The MS² fragmentation pattern under High energy Collision Dissolution (HCD) was acquired using the optimized normalized collision energy (NCE). The multistage fragmentation pathways MS_n (n= 2,3) were acquired by selecting ions of interest under Collision Induced Dissolution (CID). The fragmentation patterns obtained from both modes were compared to a known standard, diethylaminopretadalafil. The accurate mass measurement of each product ion was determined with the aid of Mass Frontier 7.0 software. Mass error of all product ion has been limited to within \pm 5.0 ppm tolerance.

Results

The fragmentation patterns for both unknown compound and diethylaminopretadalafil were compared. Under MS² experiment in HCD mode, the unknown compound showed an intense signal at m/z 492, which is more than diethylpropylaminopretadalafil at m/z 464 by 28 Dalton, indicating the presence of $-C_2H_4-$ to the known standard structure. Similarly, the difference has been observed for the fragment ion at m/z 393 of the unknown compound and the fragment ion at m/z 379 of diethylaminopretadalafil. The common ion at m/z 334 is a result of the cleavage of C-1, C-3 and nitrogen [2]. Under CID mode, the MS² spectrum of the unknown compounds showed fragmentations at m/z 393 and 292, both were 28 Dalton higher than the fragment ions at m/z 379 and 278 derived from diethylaminopretadalafil. The most intense ions at m/z 393 and 379 for both compounds can be formed from the cleavage of amido bridge via N-C bond and followed by the new N-C bond formation with a propyl and ethyl respectively. This is a good indication to rule out the possibility of having asymmetrical dialkyl groups linked to terminal amine. Two common ions have been observed at m/z 336 and 264. The first one can due to the loss of the labile *N*-propyl and *N*-ethyl for the unknown compound and diethylaminopretadalafil respectively. Meanwhile, the smaller ion at m/z 264, consisting of indole- benzodioxole substructure can be formed from the loss of the ethylmethyl ester-*N*-alkyl linked moiety.

Under MS³ experiment, the precursor ion at m/z 393 produced three major ions at m/z 336, 292 and 264. Similar fragmentation pattern had also been observed for ion at m/z 279 derived from diethylaminopretadalafil to yield fragments ions at m/z 336, 278 and 264. The difference of 12 Dalton between ions at m/z 292 and 278 supports the formation of N-C(alkyl) bond as explained above. The common fragmentation pathways for ions at m/z 336 and 264 had been observed in both MS³ experiments. The former produced an intense signal at m/z 130 due to cleavage of the ethyl methyl ester and benzodioxole from the precursor, forming the stable methylated indole. Lastly, the precursor ion at m/z 264 gave rise to ion at m/z 234 as a result of the removal of HCHO from the benzodioxole group. Subsequent neutral loss of one CO resulted in the ion at m/z 206.

Conclusions

Dipropylaminopretadalafil has been characterized by a series of comparison studies with the known diethylaminopretadalafil under High energy Collision Dissolution (HCD) and Collision Induced Dissolution (CID). The similarities and differences observed by these two compounds have been explained through their fragmentation pathways with the aids of reference ions.

Novel Aspect

In the absence of NMR facility and limited sample amount, the application of HRMS with MS_n capability is able to provide important information on the fragmentation pathways through comparison with a suitable known reference compound.

References

1. Bakota EL, Kelly AT, Walterscheid JP, Phatak DR, J. Anal. Toxicol. 41 250-255 (2017)
2. Huang YC, Lee HC, Lin YL, Tsai CF, Cheng HF, Food Addit. & Contam. Part A33 953-958 (2016)

WP-70 / EVALUATION OF SMALL MOLECULE DRUGS RETENTION IN FLUORINATED SILICA NANOPARTICLE STABILIZED PICKERING EMULSION USING ESI_MINE™ PLATFORM

Xin Li⁽¹⁾ - *Jo Francis*⁽²⁾ - *Colin Pilkington*⁽¹⁾ - *Alex Clay*⁽¹⁾ - *Clive A. Smith*⁽¹⁾

Sphere Fluidics Limited, R&D, Cambridge, United Kingdom (1) - GlaxoSmithKline (GSK), R&D, Stevenage, United Kingdom (2)

Keywords: ESI-Mine™, Pickering emulsion, Picodroplet, Microfluidics

Introduction

Microfluidic Water-in-Oil (W/O) picodroplet technology has been emerging as a useful tool for single cell and biological macromolecule characterizations. However, poor retention in polymeric surfactant stabilized picodroplets has limited its applications in small molecules related analyses.

Methods

ESI-Mine™ is a novel, high-throughput, mass spectrometry (MS) based platform developed for analyzing the contents of microfluidic picodroplets. Fluorinated silica nanoparticle (f-SiNP) stabilized W/O Pickering emulsion is reported capable of retaining a hydrophobic fluorescent dye (resorufin) in picodroplets. Here, we profiled the retention potency of non-fluorescent small molecules with a logP range from -2 to 4 in Pickering picodroplets using ESI-Mine™.

Results

Pickering emulsion picodroplets (600-700 pL) were generated on a flow focusing 60X60 Pico-Gen™ biochip. Af-SiNP(50 nm diameter) suspension in Novec™ 7500 (5%, w/w), prepared following Stöber's [1] and Pan's approaches [2], was used as the continuous phase, and an aqueous solution with/without 10 μM selected small molecule as the breaking phase. 'Occupied' picodroplets contain the selected small molecule, but 'Empty' do not. Picodroplet mixture of 'Occupied' and 'Empty' were reinjected onto a MS emitter microfluidic biochip interfacing with a Perkin Elmer Axlon 2 mass spectrometer, *i.e.* ESI-Mine™. Total ion current (TIC) signal of the picodroplets were recorded as a series of peaks. Extracted ion current (EIC) signal only shows 'Occupied' picodroplet peaks if the selected small molecule is retained inside 'Occupied' picodroplets, but in both 'Occupied' and 'Empty' if the selected small molecule cannot be retained. Eight small molecules of a logP range from -2 to 4 were selected, and analyzed on ESI-Mine™ platform.

Conclusions

We employed and demonstrated ESI-Mine™ platform as a useful tool to evaluate the retention potency of

non-fluorescent small molecules inside f-SiNP stabilized Pickering emulsion picodroplets. Improved retention potency of small molecules inside f-SiNP stabilized Pickering emulsion picodroplets was observed compared to conventional polymeric surfactant stabilized picodroplets.

Novel Aspect

Novel aspects include MS analysis of Pickering emulsion picodroplets, and evaluating retention potency of non-fluorescent small molecules inside Pickering picodroplets.

References

1. Stöber W., Fink A., J. Colloid Interface Sci., 26, 62-69 (1968).
2. Pan M., Rosenfeld L., Kim M., Xu M., Lin E., Derda R., Tang S. K. Y., ACS Appl. Mater. Interfaces, 6, 21446–21453 (2014).

WP-71 / GENERIC METHOD FOR PARALLEL ARTIFICIAL LIQUID MEMBRANE EXTRACTION OF NONPOLAR BASIC DRUGS FROM HUMAN PLASMA

Marthe Lid - Astrid Gjelstad - Stig Pedersen-Bjergaard

Pharmaceutical Institute, University of Oslo, Oslo, Norway

Keywords: Drug analysis, liquid phase microextraction (LPME), sample preparation, automation, green chemistry.

Introduction

Through parallel artificial liquid membrane extraction (PALME) it is possible to perform miniaturized liquid-liquid extraction in a 96-well plate system without the need for solvent evaporation and sample reconstitution prior to LC-MS [1-3]. Thus, with PALME, liquid-liquid extraction can be downscaled, simplified and automated. The aim of this project was to propose, develop, and test a generic method for PALME of monobasic drugs from human blood plasma with $\log P > 2$.

Methods

The setup consists of a commercially available 96-well sample plate and filter plate. Samples of human blood plasma were mixed with carbonate buffer and methanol (total volume 250 μl) and were pipetted into wells in the sample plate. Liquid membranes of 3 μl organic solvent (1 % trioctylamine in dodecylacetate) were established by pipetting the mixture on filters in the filter plate. Acceptor solutions of 50 μl aqueous formic acid (20 mM) were pipetted in the filter plate. The two plates were clamped together and agitated for 45 minutes at 900 rpm. Acceptor solutions were analyzed by UHPLC-MS/MS.

Results

Data on extraction of 46 analytes have shown recoveries $>40\%$ for most analytes. Preliminary work on linearity, precision and accuracy on a selected number of analytes complied with the requirements of international guidelines for bioanalytical method validation, and a full validation is currently in progress.

Conclusions

A generic PALME method was proposed for monobasic drugs with $\log P > 2$, and tested for extraction from human blood plasma. The preliminary results show potential and complies with international guidelines for bioanalytical method validation. The generic PALME method provided a simple workflow (with potential for automation), efficient sample cleanup for LC-MS analysis [1], and reduced the consumption of organic solvent to 3 μl per sample.

Novel Aspect

This study is the first attempt to develop generic methods for parallel artificial liquid membrane extraction

(miniaturized liquid-liquid extraction in 96-well format).

References

1. Ask, K. S., Bardakci, T., Parmer, M. P., Halvorsen, T. G., Øiestad, E. L., Pedersen-Bjergaard, S., Gjelstad, A. 2016 *Journal of Pharmaceutical and Biomedical Analysis* 129, pp. 229-236
2. Gjelstad, A., Rasmussen, K. E., Parmer, M. P., Pedersen-Bjergaard, S. 2013 *Bioanalysis* 5, pp. 1377-1385
3. Vårdal, L., Askildsen, H. -M., Gjelstad, A., Øiestad, E. L., Edvardsen, H. M. E., Pedersen-Bjergaard, S. 2017 *Journal of Chromatography B: Analytical Technologies in the Biomedical and Life Sciences* 1048, pp. 77-84

WP-72 / APPLICATION OF LINEAR EQUATIONS DECONVOLUTION ANALYSIS (LEDA) ALGORITHM IN ION TRAP MASS SPECTROMETRY: RESOLUTION OF ISOMERIC PAIRS OF CARBONIC ANHYDRASE AND CYCLOOXYGENASE HYBRIDS INHIBITORS

Donato Squillaci - Marta Menicatti - Silvia Bua - Claudiu T. Supuran - Gianluca Bartolucci

University of Florence, NEUROFARBA Dept., Sezione di Scienze Farmaceutiche, Sesto Fiorentino (FI), Italy

Keywords: linear equations of deconvolution analysis, LEDA, energy resolved mass spectrometry, isomers resolution, ion trap

Introduction

The separation of a mixtures of isomers by chromatography generally requires lot of time and it is isomers pair-specific. To overcome this issue, LC-MS/MS methods, together with a mathematical algorithm that is able to identify the isomers present in the sample without their chromatographic separation, were developed [1-2]. This approach was applied to MS/MS acquisitions obtained with both triple quadrupole and ion trap mass analyzers.

Methods

The LC-MS/MS analysis were performed using both triple quadrupole (QqQ) and ion trap (IT) systems, either equipped with ESI source, operating in positive-ion mode. The chromatographic parameters employed to analyze the samples were finely tuned to minimize the run time. Moreover, a matrix of linear equations (LEDA), based on abundance ratios of product vs precursor ions, was used. This approach was verified to 6 isomers pairs of NSAID-coumarin hybrid inhibitors [3].

Results

The study was carried out on standard solutions of analytes and their mixtures in the concentration range of 10 to 100 ng mL⁻¹. Using the reported conditions, a co-elution of the studied isomers were observed. Therefore, in case of a isomers mixture, LEDA must be able to resolve the co-eluting peak determining the relative concentration of each isomer eventually present. However, some issues were occur in the QqQ results, that showed in some cases lack of precision. These problems were mainly due to the difficult evaluation of the abundance of the precursor ion before its decomposition. To avoid the early decay of the precursor ion, the MS/MS experiments in IT were carried out. The fragmentation mechanism in IT allows a better management of the ions (precursor/products) ensuring their correct abundance rating. Indeed, the different energy transfer during collision process, lead to the formation of product ions directly by precursor ion. However, this characteristic often produce a poor number of product ions, which contrasts the LEDA application.

Conclusions

The quali-quantitative data obtained from both tested MS analyzer demonstrated the LEDA reliability in the isomers speciation.

Novel Aspect

The LEDA approach has the advantage that isomers can be quantified without the need of LC separation or additional specialized ion mobility instrumentation.

References

1. Menicatti M., Guandalini L., Dei S., Floriddia E., Teodori E., Traldi P., Bartolucci G., Rapid Commun. Mass Spectrom., 30, 423-432 (2016).
2. Menicatti M., Guandalini L., Dei S., Floriddia E., Teodori E., Traldi P., Bartolucci G., Eur. J. Mass Spectrom., 22, 235-243 (2016).
3. Bua S., Di Cesare Mannelli L., Vullo D., Ghelardini C., Bartolucci G., Scozzafava A., Supuran C. T., Carta F., J. Med. Chem., 60, 1159-1170 (2017)

WPS-S04 - Food integrity, authenticity & traceability

Chairs: Ruud Peters, Michele Suman

WP-73 / PIGMENT PROFILING IN OLIVE OILS BY ULTRA-HIGH PERFORMANCE LIQUID CHROMATOGRAPHY – ATMOSPHERIC PRESSURE CHEMICAL IONIZATION-HIGH RESOLUTION MASS SPECTROMETRY

Ane Arrizabalaga-Larrañaga⁽¹⁾ - *Pilar Rodriguez*⁽²⁾ - *Mireia Medina*⁽²⁾ - *Francisco Javier Santos*⁽¹⁾ - *Encarnación Moyano*⁽¹⁾

University of Barcelona, Chemical Engineering and Analytical Chemistry, Barcelona (1) - Government of Catalonia, Agri-Food Laboratory, Barcelona (2)

Keywords: edible oils, adulteration, food fraud, pigments, dyes

Introduction

The color is considered an important characteristic for consumers to evaluate olive oils quality. Pigments such as carotenoids and chlorophylls are the main responsible of this property but the color changes from green to yellow-brown due to the production processes or inadequate storage. In order to regreen the olive oils, instead of adding 20-30% of extra virgin olive oil, as it is established, different dyes or oil mixtures are used to adulterate [1,2].

Methods

A UHPLC-APCI-HMRS (Q-Orbitrap) method was developed to quantify chlorophylls, chlorophylls derivatives and carotenoid pigments and to detect adulterated olive oil samples due to the addition of dyes. A SPE procedure was performed using Silica SPE cartridges to clean-up the sample extracts. The UHPLC-APCI-HRMS method was achieved using an Accucore C18 column and working in positive ion mode and full scan acquisition mode (range 200-1000 m/z , Rs 70,000).

Results

The chromatographic separation of pigments was achieved in less than 15 min under gradient elution with a quaternary mobile phase. All pigments were ionized with APCI in positive ion mode. Chlorophyll and its derivatives generated the protonated molecule $[M+H]^+$ as base peak of the mass spectra with no significant in-source CID fragment ions. Concerning carotenoids, most of them also generated $[M+H]^+$ as base peak with significant in-source CID fragment although, in the case of lutein, the in-source CID fragment ion $[M+H-H_2O]^+$ dominated the mass spectrum. Moreover, tandem mass spectrometry has been used to obtain structural information for further confirmation of these analytes in olive oil samples.

The UHPLC-APCI-HRMS/MS method was applied achieving a pigment profile (relative concentrations

between the different pigments) useful for the detection of adulterated olive oils. Additionally, the pigment profile of different potential fraudulent colorants was also studied in order to identify further adulteration evidences.

Conclusions

The proposed UHPLC-APCI-HRMS/MS method showed good quality parameters and allowed the accurate pigment profile of olive oil samples. Moreover, this method provided enough evidences for the detection of adulterated olive oil by the addition of dyes.

Novel Aspect

This is the first time that UHPLC-APCI-HRMS is applied for the olive oil pigment profiling to detect a fraud due to addition of dyes.

References

1. Fang M., Tsai C., Wu G., Tseng S., Cheng H., Kuo C., Hsu C., Kao Y., Shih D. Chiang Y., Food additives and contaminants, 8, 157 (2015)
2. Roca M., Gallardo-Guerrero L., Minguez-Mosquera I., Gandul B., Journal of Agricultural and food chemistry, 58, 51 (2010)

WP-74 / ISOTOPIC MAPPING OF MILK AND GROUNDWATER IN THE PRODUCTION AREA OF PARMIGIANO REGGIANO FOR FOOD TRACEABILITY

Marta Boito⁽¹⁾ - ***Paola Iacumin***⁽¹⁾ - ***Nives Ogrinc***⁽²⁾

University of Parma, Department of Chemistry, Life Sciences and Environmental Sustainability, Parma, Italy
(1) - Jožef Stefan Institute, Ljubljana, Slovenia (2)

Keywords: isotopes, milk, origin, Parmigiano Reggiano, feeding regime

Introduction

Determination of the geographical origin of foodstuffs is becoming of increasing interest to consumers and producers, since it may be used as a criterion for certifying quality and authenticity.

The Protected Designation of Origin (PDO) trademark has been assigned to numerous local products based on their area of origin. In order to obtain this designation, the raw materials must have been produced and processed in the specific region from which the product gets its name. [1][2]

Methods

To make a check at the base of the production of "Parmigiano Reggiano" we began to collect monthly samples of groundwater, milk, hay and fodder from ten cattle-sheds in the production area. Sampling will last for at least a year. The stable isotope composition of the elements, such as O, H, C, N, have been evaluated with the aim of characterizing the milk of the Consortium areas and defining the isotopic fractionation of water ingested/milk produced and milk produced/food supplied.

Results and Conclusions

From the first results obtained it was possible to notice a differentiation of the isotopic values of O and H between plain and mountain milk. Even, the hay produced locally by breeders is different compared to the one purchased by import. Even, these greater characterizations could increase the added value of one producer's cheese.

Novel Aspect

The results of the present study could be helpful for consumers, the food industry and government

regulatory agencies as it can prevent fraudulent labelling of organic food. [3]

References

1. D. Sacco, M. A. Brescia, A. Sgaramella, G. Casiello, A. Buccolieri, N. Ogrinc, A. Sacco, Discrimination between Southern Italy and foreign milk samples using spectroscopic and analytical data, *Food Chemistry* 114 (2009) 1559–1563
2. Dotsika E., Bulleri E., Brunella R., Nudda A., Mele M., Battaglini R., Poutoukis D. and Chantzi P., Fingerprinting the origin of sheep milk products based on H-O-C-N isotopes, *Integrative Food, Nutrition and Metabolism*, 2017
3. Ill-Min Chung, Inmyoung Park, Jae-Yeon Yoon, Ye-Seul Yang, Seung-Hyun Kim, Determination of organic milk authenticity using carbon and nitrogen natural isotopes, *Food Chemistry* 160 (2014) 214-218

WP-75 / WINE TRACEABILITY: "VIGNETO ITALIA" PROJECT

Claudio Brilli - Benedetta Bronzi

Ruffino srl, Laboratory, Pontassieve (Florence), Italy

Keywords: wine, traceability, ICP-MS, Chianti

Introduction

The project aims to classify wines on the basis of different geographical origin of the grapes. 246 and 393 grapes from known provenance are collected during 2013 and 2014 harvests. Then they are turned into wine by a vinification procedure that simulates the industrial process. Elemental composition of these samples is analyzed with ICP-MS and the resulting data are used to train a classification model that discriminates Chianti from other wines.

Methods

Two quantitative procedures are used depending on the concentration level of the analytes. Most elements (10 ppb - 1000 ppm) are analyzed directly after 1:20 dilution and standard nebulization into the ionization chamber. For ultratrace elements (conc. < 100 ppt) the sample is diluted 1:5 after mineralization of the matrix; in this case the sample introduction system (Apex-Q with Spiro TMD) is designed to increase analytical sensitivity and to reduce interferences.

Results

After statistical pre-processing of the data, three classification models are obtained through chemometric elaboration of ICP-MS analysis results: Linear Model and Similarity Model are based on supervised pattern recognition methods; Non Linear Model is based on Artificial Neural Networks (ANNs) which is one the most promising techniques to cope with the problem of automatic and intelligent classification. An additional model (Consensus) calculates the average between the results of the previous three models. The predictive models are validated by testing 134 wine samples with known geographical origin: the Consensus Model shows a non-error rate (NER%) of 89,9% with 18,7% of non-classified samples. Wine adulteration is also evaluated by analysing Chianti and non-Chianti samples blended at different rates: the results demonstrate that Consensus Model is able to detect down to (at least) 25% addition of non-Chianti wine to a Chianti sample, since it classifies such mixture as non-Chianti.

Conclusions

The predictive model Consensus shows a satisfying percentage of correct classification. Some regions (such as Sicily) are difficult to differentiate from the Chianti area. But many regions with high wine production (Sardinia, Piedmont, Puglia, Marche, Campania, Basilicata, Umbria) are correctly identified. Besides the total overlapping of samples from vintages 2013 and 2014 suggests that correlation is not

affected by vintage.

WP-76 / UHPLC-ESI-MS/MS FOR THE DETERMINATION OF POLYPHENOLS IN SPANISH PAPRIKA (CAPSICUM ANNUUM L.). APPLICATION TO PROTECTED DESIGNATION OF ORIGIN (PDO) CLASSIFICATION BY CHEMOMETRICS.

Guillem Campmajó - Sergio Barbosa - Oscar Núñez - Javier Saurina - Lluís Puignou

University of Barcelona, Department of Chemical Engineering and Analytical Chemistry, Barcelona

Keywords: Polyphenols; Spanish Paprika; Food Authentication; UHPLC; Tandem Mass Spectrometry

Introduction

Polyphenols are a large family of aromatic secondary metabolites well-known for their health beneficial effects and for contributing to sensorial attributes [1]. These compounds are mainly distributed in food of plant origin such as paprika, which is a dried and ground spice obtained from red pepper. In Spain, only La Vera and Murcia varieties are distinguished with protected designation of origin (PDO) and their polyphenolic content may allow to classify them.

Methods

A UHPLC-ESI-MS/MS (QqQ) method was developed to determine 37 polyphenols in paprika samples. A simple sample extraction procedure using H₂O:ACN (20:80, v/v) was employed. The chromatographic separation was optimized using an Ascentis Express C18 reversed-phase (10 cm × 2.1 mm, 2.7 μm) fused-core column (analysis time of 30 min). Paprika polyphenolic contents were subjected to principal component analysis (PCA) using PLS_Toolbox 7.8.2 (Eigenvector Research).

Results

The developed UHPLC-ESI-MS/MS method showed satisfactory limits of detection (down to 0.012 μg/L), linearity ($r^2 \geq 0.995$), run-to-run and day-to-day precisions ($RSD \leq 20\%$), and trueness (relative errors below 15%). 111 paprika samples from La Vera and Murcia (Spain), and from Czech Republic were analyzed, and the polyphenolic content quantified by external calibration. Chlorogenic, ferulic and p-coumaric acids, hesperidin, homoplantagin, rutin, and vanillin were found in all the samples. Polyphenolic composition was considered as a source of potential descriptors to be exploited for the classification of paprika according to the PDO and region of production. The plot of scores showed a successful separation in accordance with both PDO and production region. Moreover, a discrimination between Murcia and Czech Republic paprika varieties (sweet and spicy) was also possible. The PCA loading plots revealed that homoplantagin, betulinic acid and kaempferol were among the most relevant polyphenols to achieve this discrimination.

Conclusions

A simple and reliable UHPLC-ESI-MS/MS method for the quantitation of 37 polyphenols in paprika samples was developed. The characterization and classification of paprika samples according to the PDO and production region was achieved by employing polyphenolic content as chemical descriptors for PCA. This method will be useful to guarantee PDO attributes and for the prevention of frauds in paprika samples.

Novel Aspect

Fast UHPLC-ESI-MS/MS method for the determination of polyphenols and the classification of paprika samples in the authentication of PDOs.

References

1. Lucci P., Saurina J., Núñez O., TrAC Trends Anal. Chem., 88, 1 (2017).

WP-77 / A NOVEL MULTI-PLATFORM HIGH RESOLUTION MASS SPECTROMETRY NON-TARGETED APPROACH FACING EXTRA VIRGIN OLIVE OIL ADULTERATION

Daniele Cavanna^(1,2) - ***Kamila Hurkova***⁽³⁾ - ***Zbynek Dzuman***⁽³⁾ - ***Andrea Serani***⁽⁴⁾ - ***Chiara Dall'Asta***⁽²⁾ - ***Jana Hajslova***⁽³⁾ - ***Michele Suman***⁽¹⁾

Advanced Laboratory Research, Barilla G. e R. Fratelli S. p. A., Via Mantova 166-43122 Parma (Italy) (1) - Department of Food and Drug, University of Parma, Parco Area delle Scienze 95/A-43124 Parma,(Italy) (2) - University of Chemistry and Technology, Prague, Department of Food Analysis and Nutrition, Technicka 3, 166 28 Prague 6 (Czech Republic) (3) - COTECA Srl Consulenze Tecniche agroindustriali - 56121 Pisa (Italy) (4)

Keywords: Non-targeted mass spectrometry; Extra Virgin Olive Oil; Food frauds; LC-HRMS

Introduction

High resolution mass spectrometry was used in the past for the detection of frauds on extra virgin olive oil products [1] [2] [3]. However, at the moment, scientific studies are not largely focused on the detection of soft refinements of Olive Oils, even if the dilution of EVOO with these samples is probably one of the most common frauds applied on the market.

This inter-laboratory work would like to identify some chemical markers responsible of this process.

Methods

Refined oils (deodorized and deacidified) were created on a laboratory scale starting from low quality olive oils and were analyzed together with a set of pure EVOO samples.

Sample preparation was executed according to previous works [4] [5].

UHPLC- HRMS analyses were performed with a SCIEX® triple TOF 6600 analyzer and with a Thermo® Q-Exactive Orbitrap analyzer.

Instrumental conditions were the same for all the instruments involved in the project.

Results

Raw data were processed with PeakView® (for triple TOF data) and Compound Discoverer® (for Orbitrap data); chemometric data analyses were performed with Markerview® (Sciex) and SIMCA® software (Umetrics).

Different PCA models were created and the scores plots clearly highlighted a separation between the pure Extra Virgin Olive Oil samples and the soft refined samples.

The features responsible of these clusterizations were selected and an attempt of compound identification was performed.

A group of molecules related to the adulteration processes were identified: they were at least upregulated in the refined oils or even only detected in these fraudulent samples and in their mixtures.

The markers selected applying the same experimental design with different instruments in different laboratories were compared.

Conclusions

This inter-laboratory research represents a robust untargeted metabolomic approach able to identify specific compounds that can be considered markers of soft refinement processes in Extra Virgin Olive Oils.

The detection of these molecules (performed with easier “target” methods) will be a proof of fraudulent issues in commercial samples.

Novel Aspect

Chemical markers for this important fraud in EVOO have not been largely studied yet. The multi-platform non-targeted approach could also be used for retrospective analyses.

References

1. Gil-Solsona R., Raro M., Sales C., Lacalle L., Díaz R., Ibáñez M., Beltran J., Vicente Sancho J., Hernández F. J., *Food Control*, 70, 350-359 (2016)
2. Kalogiouri N. P., Alygizakis N. A., Aalizadeh R., Thomaidis N. S., *Analytical and Bioanalytical Chemistry*, 408, 7955-7970 (2016)
3. Di Girolamo F., Masotti A., Lante I., Scapatucci M., Calvano C. D., Zambonin C., Muraca M., Putignani L., *International Journal of Molecular Sciences*, 16, 20896-20912 (2015)
4. Vaclavik L., Cajka T., Hrbek V., Hajslova J., *Analytica Chimica Acta*, 645, 56-63 (2009)
5. Hurkova K., Rubert J., Stranska-Zachariasova M., Hajslova J., *Food Analytical Methods*, 10, 1317-1327 (2017)

WP-78 / SHRIMP FRAUD: EVALUATING CHARACTERISTIC METABOLOMIC FEATURES FOR AUTHENTICATION BY HIGH AND UNIT MASS RESOLUTION SPECTROMETRY

Olivier Chevallier⁽¹⁾ - Niladri Chatterjee⁽²⁾ - Ewa Wielogorska⁽³⁾ - Connor Black⁽⁴⁾ - Christopher Elliott⁽⁴⁾

Queens University Belfast, Mass Spectrometry Core Unit, Belfast (1) - ICAR, Central Institute of Fisheries Technology, Cochin (2) - Queens University Belfast, School of Pharmacy, Belfast (3) - Queens University Belfast, Institute for Global Food Security, Belfast (4)

Keywords: Food fraud, Food authenticity, LC-HRMS, LC-MS/MS, chemometric

Introduction

Unfortunately fishery products occupy the second highest ranking position among commodities that are at most risk of food fraud. Among these, marine shrimps and prawns accounts for more than 17% of the global seafood consumption. As per Article 35 of the EC regulation 1379/2013, commercial designation of the species, production method and geographical origin have to be provided to consumers.

Methods

Authentic samples were collected either directly from aquaculture farms or through local supermarket chains. Six species of shrimps namely Tiger Prawn, King Prawn, Indian Prawn, Pink Speckled Shrimp, Argentinian Red Shrimp and Red Prawn were included in the study. Among Tiger Prawns, wild caught prawns from India and Madagascar; and farmed prawns from Vietnam were sources as were farmed King Prawns were from India, Thailand, Vietnam and Honduras.

Results

Methanolic extracts of freeze-dried samples were analysed in full scan positive and negative ionisation mode using a Waters Acquity I-Class UPLC coupled to a Xevo G2S QToF. Chromatographic alignment, peak picking and feature identification of the acquired data were carried out on Progenesis Q1 software followed by multivariate analysis in SIMCA-P software. OPLS-DA chemometric models could discriminate between species, geographical origin and production method. The R² and Q² values for all the models were >0.9. Further, characteristic metabolomic features for each species were determined by group wise comparison of the data. These features were evaluated in a targeted, triple quadrupole LC-MS/MS method and at least one exclusive biomarker for each species was selected. Subsequently, a confirmatory LC-MS/MS method for authentication of species was developed. An OPLS-DA model built using the characteristic features of Tiger Prawns and King Prawns was shown to discriminate between the geographical origins of these species.

Conclusions

The data generated gives a strong indication that the routine testing for species authenticity as well as the geographical origin of the prawns can be achieved using triple quadrupole LC-MS/MS.

Novel Aspect

This innovative methodology based on LC-HRMS and LC-MS/MS can present an alternative to DNS testing or isotope based assays

WP-79 / OPTIMIZATION, VALIDATION AND APPLICATION OF GC-FID AND GC-IRMS METHODS FOR THE DETERMINATION OF FATTY ACIDS WITHOUT DERIVATISATION

Andrija Ciric - Jasmina Masten - Lidija Strojnik - Nives Ogrinc

Jozef Stefan Institute, Department of Environmental Science, Ljubljana, Slovenia

Keywords: Fatty acids, FFAP, stable isotope ratio

Introduction

The fatty acid (FA) composition of food is a recurring point of attention in government policy concerning nutrition and health, since they are source of energy, components of biological membranes, precursor for many different molecules and as transporters for vitamins such as A, D, E and K [1]. Further the stable isotope composition of carbon and hydrogen in fatty acids can be used to assess the geographical origin and authenticity of food products [2]. The aim of this presentation was to develop a simple, robust, and derivatization-free GC method suitable for the routine analysis of the major FA components and their isotope compositions.

Methods

All analyses were performed using an Agilent model Hewlett-Packard 6890 gas chromatograph equipped with a flame ionization detector (GC-FID), while for stable isotope analysis an isotope ratio mass spectrometer IsoPrime (GV Instrument) with a combustion interface (GC-C-IRMS) was used. Separation was achieved using a Phenomex ZEBRON FFAP (30 m × 0.25 mm with film thickness of 0.25 μm) capillary column. The temperature program was as follows: 130°C to 255°C at rate of 20°C/min. Helium was used for the carrier gas at a flow rate of 2.6 mL/min.

Results

The separation of the FA was optimized by adjusting both the initial oven temperature and the temperature ramp. It was noticed that a faster temperature rate resulted in better resolution between palmitic and palmitoleic acid. The linearity of the method was verified between 50 - 200% of the injection volume as prescribed by the method. The results show an acceptable RSD% and mass deviation for all 16 major and minor constituents and confirm the linearity of the method. Method precision was checked in terms of repeatability and reproducibility. When the method was optimized using GC-FID, it was also applied to determine the stable isotope composition of C by GC-C-IRMS.

The stable carbon isotopic values are expressed in the usual notation as the per mill (‰) deviation of the isotope ratio of a sample relative to that of a standard. The $\delta^{13}\text{C}$ values of the 16 main FA vary from -30.9 to -28.6 ‰. The developed method was further tested on milk and oil samples.

Conclusions

The method for determining the FA profile of different food samples has been optimized for the nitroterephthalic acid phase. results showed that the method has acceptable selectivity, linearity precision, detection and quantification limits, robustness and stability. The method is suitable for the routine determination of the amount and isotope composition of FA in vegetables and food of animal origin with a distinct FA profile.

Novel Aspect

Determining the fatty acid profile of different food samples without the need for derivatisation. The method can be used for food nutritional labelling, authentication, traceability and quality control.

References

1. Serhan C. N., Chiang N., Van Dyke T. E., Nat Rev Immunol., 8, 349-361 (2008).
2. Liu X., Xu S. P., Wang J. H., Yuan J. P., Guo L. X., Li X., Huang X. N., Anal. Bioanal. Chem., 388, 723-731 (2007).

Acknowledgment

The study is part of two EU projects: MASSTWIN (H2020, grant agreement no. 692241) and ERA Chair ISO-FOOD (FP7, grant agreement no. 621329).

WP-80 / MULTI-CLASS LC-MS AND GC-MS METHODOLOGIES TO DISCRIMINATE EXTRA VIRGIN OLIVE OILS FROM DIFFERENT PROTECTED DESIGNATIONS OF ORIGIN

Lucía Olmo-García⁽¹⁾ - *Karin Wendt*⁽²⁾ - *Aadil Bajoub*⁽³⁾ - *Carmen María Sánchez-Arévalo*⁽¹⁾ - *Arsenio Muñoz de la Peña*⁽⁴⁾ - *Carsten Baessmann*⁽²⁾ - *Alegría Carrasco-Pancorbo*⁽¹⁾ - *Alberto Fernández-Gutiérrez*⁽¹⁾

University of Granada, Department of Analytical Chemistry, Granada (1) - Bruker Daltonik GmbH, Bruker Daltonik GmbH, Bremen (2) - National School of Agriculture, Department of Basic Sciences, Meknès (3) - University of Extremadura, Department of Analytical Chemistry, Badajoz (4)

Keywords: virgin olive oil; protected designations of origin; LC-MS, GC-MS; chemometrics.

Introduction

Product labeling has gained considerable attention over the last years, as a means to both provide product-specific information and reduce the quality uncertainty faced by consumers. Implementation of Protected Designations of Origin (PDOs) is one of the most prominent differentiation strategies used in olive oil market. They are often perceived as valuable tools that promote specific attributes of the oil linked to its geographical provenance.

Methods

126 oil samples from 6 Mediterranean PDOs were analyzed by LC-MS and GC-MS combined to chemometrics. The extracts were eluted in LC using a C18 (2.1x 100 mm, 1.8 µm) column, with acidified water and acetonitrile and a flow gradient (0.4-0.6 mL/min) at 40°C. The derivatized extracts were injected in GC, using a BR-5 column with a T gradient from 150 to 320°C (4°C/min rate). Both systems were coupled to a Compact™ QTOF MS (Bruker) by an ESI interface for LC and an APCI source for GC.

Results

Samples from six different Mediterranean PDOs (Meknès and Ouazzane (Morocco), Priego de Córdoba and Baena (Spain), Kalamata (Greece) and Toscana (Italy)) were collected and analyzed. Sample prep. (unselective 3-steps liquid-liquid extraction) and experimental conditions facilitated the determination of a great number of minor compounds (phenolic compounds, triterpenic acids and dialcohols, sterols, tocopherols, free fatty acids, etc.) within a single run.

Non-targeted and targeted approaches were used to offer maximum coverage of the olive oil metabolome's chemical space in a 1st step, and the possible validation of the identified markers afterwards. Data treatment (PCA, PLS-DA) -done by using MetaboScape®- led to a noticeable discrimination among the six evaluated PDOs taking into account the data coming from both platforms. Several compounds such as elenolic acid, luteolin, oleuropein and ligstroside aglycones, and some other tentatively identified substances, were pointed out as possible PDOs distinctive markers.

Conclusions

- Multi-class LC-MS and GC-MS methodologies together with chemometrics made possible the discrimination among different PDOs, identifying potential origin markers.
- The combined use of non-targeted and targeted approaches can enhance or reinforce the outcomes of any

study.

-GC-APCI-Q TOF can preserve the pseudo-molecular ion information, which is a great advantage over the “classical” GC systems and facilitates the identification.

Novel Aspect

Combination of non-targeted and targeted approaches implying the use of powerful platforms (LC-ESI-Q TOF/GC-APCI-Q TOF) to identify potential “PDOs’ markers”.

WP-81 / DETERMINATION OF MELATONIN AND RELATED COMPOUNDS IN BEVERAGES USING LIQUID CHROMATOGRAPHY COUPLED TO LOW AND HIGH RESOLUTION MASS SPECTROMETRY

María Luisa Jiménez Medina - Lorena Manzano Sánchez - Jesús Marín Sáez - Roberto Romero González - Antonia Garrido Frenich

University of Almería, Dept. Chemistry and Physics, Almería

Keywords: Melatonin, related compounds, beverages, triple quadrupole, high resolution mass spectrometry

Introduction

Melatonin (*N*-acetyl-5-methoxytryptamine) is an endogenous indolamine hormone, synthesized from tryptophan [1] and it has a significant role in regulation of the circadian rhythm, mitigation of sleeping disorders and jet lag. Melatonin can also be found in a variety of foods including edible plants, fruits and beverages [2] in addition to isomers or related compounds, but their identification is not a straightforward process [3].

Methods

For the determination of melatonin and related compounds, a simple dilution of the beverage, including different types of beer (alcoholic and non-alcoholic), wine, aloe vera drinks and cherry juice, was performed prior determination applying liquid chromatography (LC) coupled to triple quadrupole tandem mass spectrometry (QqQ-MS/MS) or high resolution mass spectrometry analyzer as Exactive-Orbitrap.

Results

First, a nontargeted approach based on LC coupled to Exactive-Orbitrap has been used for the reliable identification of melatonin related compounds in beverages. Thus, in addition to melatonin, tryptophan ethyl ester, which has the same molecular formula than melatonin, was detected. All ion fragmentation (AIF) mode was used to get MS/MS spectra and different ions were obtained for both compounds. Then, an analytical method based on LC-QqQ-MS/MS has been validated in terms of selectivity, limit of quantification (0.01 µg/L) and precision, for the determination of melatonin and tryptophan ethyl ester in beverages. It was observed that tryptophan ethyl ester was detected at higher concentrations (from 0.2 to 90 µg/L) than melatonin (<2 µg/L) in the tested samples. In addition, the concentration of tryptophan ethyl ester was higher in alcoholic beers than in non-alcoholic, indicating that the fermentation could affect the concentration of this compound.

Conclusions

Melatonin related compounds, as tryptophan ethyl ester, has been identified in several beverages using Exactive-Orbitrap analyzer and combining the information provided by AIF mode and characteristic fragmentation pattern of the molecule. Because both compounds have similar retention time and common ions, false positives of melatonin could be reported in some beverages due to tryptophan ethyl ester was detected at higher concentrations.

Novel Aspect

The combination of several MS platforms allows for the reliable identification and determination of melatonin and related compounds in different types of beverages.

References

1. Reinholds R., Pugajeva I., Radenkovs V., Rjabova J., Bartkevic V., *Journal of Chromatographic Science*, 54, 977-984 (2016).
2. Kocadagli T., Yilmaz C., Gökmen V., *Food Chemistry*, 153, 151-156 (2014).
3. Gardana C., Iriti M., Stuknytė M., De Noni I., Simonetti P., *Journal of Pineal Research*, 57, 435-441 (2014).

WP-82 / UHPLC-ESI-HRMS CHARACTERIZATION OF TURMERIC AND CURRY SAMPLES BASED ON PROFILING AND FINGERPRINTING APPROACHES

Nerea Núñez - Oscar Núñez - Sonia Sentellas - Javier Saurina

University of Barcelona, Department of Chemical Engineering and Analytical Chemistry, Barcelona

Keywords: Curcumin; Polyphenols; Turmeric; Food Authentication; UHPLC-HRMS

Introduction

Turmeric (*Curcuma longa*) is a plant related to ginger family that has been used for centuries as a remedy in traditional Asian pharmacy due to its antioxidant, anti-inflammatory and antineoplastic activities [1]. Turmeric is also appreciated worldwide as a condiment in cuisine. Curcumin is the most relevant molecule of turmeric, providing both color and pharmacological activity [1]. Related curcuminoids and polyphenols are important as well.

Methods

The UHPLC-HRMS method used an Accela UHPLC system coupled to a Q-Exactive q-Orbitrap mass spectrometer (Thermo Fischer Scientific) equipped with a heated-electrospray ionization (H-ESI) source. Compounds were separated by reversed-phase mode using an elution gradient based on 0.1 % formic acid in water (v/v) and acetonitrile. Chromatograms were acquired in the negative mode in the range 100 to 1000 *m/z*. Data were analyzed using PLSToolbox (Eigenvector Research).

Results

Here, UHPLC-HRMS was used to analyze commercial turmeric and curry samples. The resulting data consisting on polyphenolic profiles and UHPLC-HRMS fingerprints was treated chemometrically to tackle characterization and authentication issues. Quality controls (QCs) prepared from all the samples were used to assess the instrumental variability. Principal component analysis (PCA) models from both approaches showed that curry and turmeric samples were clearly segregated. QCs appeared in a compact group in the center of the model, showing the good performance and reproducibility of the proposed method and the obtained data. Principal component 1 (PC1) mainly described quantitative patterns, with samples richer in curcuminoids to the right of the plots and those with lower contents to the left. Turmeric samples were distributed according to plant varieties, *e.g.*, Erode in the top right corner, Allepey in the bottom right corner, and Madras in the top left corner.

Conclusions

UHPLC-HRMS resulted in an excellent approach for a feasible determination of curcumin and its derivatives in food samples. The proposed method revealed a great complexity in polyphenolic composition, showing a variety of components to be elucidated. It was believed that minor curcuminoids could be relevant descriptors for classification and authentication purposes using chemometric methods.

Novel Aspect

Comprehensive UHPLC-ESI-HRMS characterization of polyphenols and curcuminoids and classification of turmeric products by PCA.

References

1. Marchiani, A. ; Rozzo, C. ; Fadda, A. ; Delogu, G. ; Ruzza, P. *Curr. Med. Chem.*, 21, 204 (2014).

WP-83 / MONOTERPENES CHARACTERIZATION IN GRAPE SPIRITS USING FAST GAS CHROMATOGRAPHY TANDEM MASS SPECTROMETRY

Mauro Paolini - S. Poli - Loris Tonidandel - Sergio Moser - Roberto Larcher

Technology Transfer Centre, Fondazione Edmund Mach, via E. Mach 1, 38010 San Michele all'Adige, Italy

Keywords: Monoterpenes, Fast GC-MS/MS analysis, Grape spirits

Introduction

Monoterpenes are among the most important aromatic metabolites to be found in a free and bound form in grapes [1]. Aromatic characteristic of grape spirits can therefore be especially attributed to this class of molecules despite the presence of countless volatile compounds produced during alcoholic fermentation [2,3]. The aim of this study was to describe the terpenic profiles responsible for the varietal aroma of a selection of grape spirits produced from 10 different varieties using an innovative fast GC-MS/MS method.

Methods

20 mL of sample was diluted to 100 mL with Milli-Q water after the addition of 2-octanol as internal standard. Monoterpenes were extracted by adsorption on a SPE cartridge and eluted with dichloromethane [4,5]. The organic phase was analysed using an Agilent Intuvo 9000 GC system coupled with an Agilent 7000 Series Triple Quadrupole MS operating in dynamic MRM mode. Separation was obtained injecting 2 μ l in split mode (1:5) into a DB-Wax UI (20 m \times 0.18 mm \times 0.18 μ m).

Results

The presented method permitted to qualify and quantify 15 different monoterpenes characteristic of the grape spirits in only 15 minutes. The compounds detected were α -terpineol, β -citronellol, geranic acid, geraniol, HO-diol I, HO-trienol, linalool oxide A, linalool oxide B, linalool oxide C, linalool oxide D, linalool, nerol, rose oxide I, rose oxide II and terpinen-4-ol. All of the measured monoterpenes were found in the 39 grape spirits except for HO-diol I and HO-trienol that were detected in only 12 and 30 samples, respectively. The correlation matrix showed a positive correlation between the amount of linalool oxides, linalool, α -terpineol and HO-diol I. A positive correlation was also observed between β -citronellol, nerol, geraniol and rose oxides.

Principal component analysis (PCA) was used to display the samples in an unsupervised pattern recognition map. The PCA showed that a group of samples was characterised by a high content of linalool oxides and linalool, while another one by a high content of β -citronellol, nerol and rose oxides.

Conclusions

The proposed method made it possible to obtain a broad monoterpene profile characteristic of grape spirit products in only 15 minutes. 15 monoterpenes were qualified and quantified in 39 grape spirit samples obtained from different grape varieties. The statistical approach highlighted that there was a positive correlation between the different compounds and allowed to discriminate the most aromatic samples on the basis of their different monoterpene profile.

Novel Aspect

The method developed can be proposed for the fast characterisation of grape spirit samples based on their

different monoterpene profile.

References

1. Maicas S., Mateo J. J., Appl. Microbiol. Biotechnol., 67, 322(2005).
2. Ferreira V., LopezR., AznarM., Analysis of Taste and Aroma, Springer, Berlin,89 (2001).
3. Francisl. L., Newton J. L., Aust. J. Grape Wine Res., 11, 114(2005).
4. Gunata Y. Z., Bayonove C. L., Baumes R. L., Cordonnier R. E., J. Chromatogr., 331, 83 (1985).
5. Boido E., Lloret A., Medina K., Farina L., Carrau F., Versini G., Dellacassa E., J. Agric. FoodChem., 51, 5408 (2003).

WP-84 / ELUCIDATION OF OLIGOSACCHARIDES IN MILK REFERENCE MATERIALS BY HILIC-MS/MS AND THE NIST TANDEM MS LIBRARIES

Connie A. Remoroza - Stephen E. Stein

Mass Spectrometry Data Center, Biomolecular Measurement Division, National Institute of Standards and Technology, Gaithersburg, MD, 20899-8362, USA

Keywords: Human Milk Oligosaccharides, Fourier Transform-ITMS, HILIC-MS-HCD, Isomers, MS library

Introduction

A very low percentage of mothers around the world practice exclusive breastfeeding. The low breastfeeding rate is the driving force for the infant formula manufacturers to produce formula that simulate human milk as closely as possible. However, the information of bioactive molecules in the NIST *SRMs for human and infant formula milk is limited. This study focuses on the method of structural elucidation of a material based mass spectral library of functional oligosaccharides using NIST reference materials.

Methods

The human milk, infant formula, goat and buffalo milk underwent sample preparation [1]. The purified oligosaccharides were analyzed on Hydrophilic Liquid Interaction Chromatography-electrospray ionization in Fusion Lumos Orbitrap-based MS using HCD & Fourier transform ion trap (FT-IT) fragmentation techniques. The spectra are extracted from LC-MS raw data and processed using in-house algorithm then analyzed using different tools including NIST Tandem mass spectral library.

Results

The high mass accuracy both for precursor and product ions from FT-IT and HCD data was sufficient to enable the assignment of oligosaccharides including novel oligosaccharide structures with different adduct types in negative/positive detection mode. Core group (Gal-Glc) with different terminal sugars, Fuc-Gal-Glc eluted prior to Sia-Gal-Glc. The elution profile of DP 5 isomers, such as lacto-*N*-fucopentaose (I, II, III) and LS-tetrasaccharide (a, b, c) displayed different retention times on HILIC that enhances oligosaccharide identification. A clear distinction in the HILIC elution pattern of isomeric oligosaccharides related to size and polarity confirmed the identifications. The base peak chromatogram of infant formula differs from human milk in terms of their glycan profiles. Infant formula contains fucolactose as well as maltodextrins with degree of polymerization 2-13 whereas human milk has over 70 bioactive oligosaccharides. The spectral library of oligosaccharides in human milk demonstrate to aid in identifying glycans and other compounds present in other biological samples such Asian goat and buffalo milk, which represents a new variety of mass spectral library. Moreover, buffalo milk contains oligosaccharides close to human milk.

Conclusions

The HILIC MS/MS and NIST Tandem MS libraries facilitate the structural assignment of 70 complex

oligosaccharides with 470 spectra including newly annotated glycan structures. Infant formula was fortified with maltodextrins and fucolactoses— a concept that can be extended to other varieties of biological fluids such as detection of glycans in other milk and human plasma samples.

Novel Aspect

A comprehensive structural and tandem MS library of annotated oligosaccharides derived from the NIST milk reference materials: a robust method of detecting other compounds in milk products.

References

1. Wu, S. et al. (2010). Journal of Proteome Research 9, 4138-4151

* SRM (Standard Reference Material)

WP-85 / IKNIFE, A NOVEL AMBIENT MASS SPECTROMETRY APPROACH FOR THE REAL-TIME IDENTIFICATION OF FOOD SAMPLES

Francesca Rigano⁽¹⁾ - ***Domenica Mangraviti***⁽²⁾ - ***Luigi Mondello***⁽²⁾

Chromaleont Srl, c/o Dipartimento ChiBioFarAm, University of Messina, Messina (1) - University of Messina, Dipartimento ChiBioFarAm, Messina (2)

Keywords: Rapid evaporative ionization mass spectrometry, electrosurgical knife, food authenticity, food safety, chemometric data analysis

Introduction

Rapid Evaporative Ionization Mass Spectrometry (REIMS) is a novel ambient MS method that can be used for the real-time identification of unknown samples, without any sample pretreatment. Such an approach is based on the analysis of the vapor produced by thermal ablation. MS profile represents an univocal fingerprinting, usable for geography evaluation and authenticity assessment of foodstuffs. For this purpose, a database of authentic samples need to be created.

Methods

The REIMS method was employed for the fast characterization of different food products (pistachio, fish, cheese, olive oil). It operates using an electrosurgical knife (iknife), which creates an aerosol that is evacuated from the sample through a transfer line into the ionization source, where a heated collision surface is located for the thermal ionization process. A multivariate statistical algorithm was validated for the real-time identification process.

Results

High value products, granted with PDO (Protected Designation of Origin) indication, were analyzed and successfully differentiated, without any geographical mismatching, thus demonstrating the applicability of the new technique in the detection of food fraud. The new technology was also applied to commercially popular and genetically similar Mediterranean Sea fish to obtain fast and accurate speciation results.

Marker compounds, responsible of sample discrimination, were determined and confirmed by conventional chromatographic techniques. These results demonstrated the capability of the iknife to provide in few seconds an holistic profile of the matrix.

Moreover, inter-laboratory studies gave proof of the reproducibility of such MS approach and transferability of the statistical model. Within this context, the present study confirmed also the potentiality of the new technique to be easy-to-use by any operator, even not expert, after the building of a reliable database and a statistical model by an analytical chemist.

Conclusions

Iknife represents a powerful tool in the preservation of food security and safety. It can be also used as a

shotgun approach to achieve a comprehensive characterization of a complex sample and, since lipids are normally the most representative detected molecules, it can be successfully employed in lipidomics. Finally, the comparison with chromatography methods made marker identification very reliable, conversely from the tentative assignment of previous works.

Novel Aspect

Among ambient methods, the iknife has been developed in close combination with a dedicated statistical software, that made the identification process very straightforward.

WP-86 / DETERMINATION OF CROCIN ISOMERS AND ISOFORMS BY MEANS UHPLC-MS/MS AS QUALITY MARKER FOR SAFFRON

Rachele Rocchi - Marcello Mascini - Manuel Sergi - Dario Compagnone - Paola Pittia

University of Teramo, Faculty of Bioscience and Technology for Food, Agriculture and Environment, Teramo

Keywords: saffron, crocins, UHPLC-MS/MS, process marker, multivariate analysis

Introduction

Saffron is a precious spice known as “red gold” obtained from the dried stigmas of the *Crocus Sativus* L., a plant of the Iridaceae family. Its importance is related to the peculiar quality and sensory properties, in particular aroma, colour and taste that are mainly due to the presence of safranal, crocins and picrocrocin, respectively [1]. Recent studies have highlighted the different beneficial health values beyond its colouring and flavouring abilities [2].

Methods

Forty-two saffron samples of different origin, drying process and age were collected. An UHPLC equipment composed by a Nexera LC20AD XR apparatus (Shimadzu, Tokyo, Japan) coupled with a 4500 Qtrap mass spectrometer (Sciex, Toronto, Canada) equipped with a Turbo V ESI source, was used for analysis. The target analytes were separated using a Kinetex C18 column in 9 minutes. The results were analysed by means of principal component analysis (PCA) using XLSTAT2016.

Results

A preliminary experimental design was applied to optimize the extraction of crocins. UHPLC-MS/MS conditions were set to obtain the best analytical performances in terms of sensitivity and selectivity; the selected experimental conditions allowed to determine all crocin isomers and isoforms in a single run. The crocin content in the samples was significantly different and resulted affected by process, age and origin, with a clear separation between the mild and high thermally processed samples. In particular results indicate that, increasing the drying temperature, a degradation of the crocin isomers with high number of glucose moieties occurs. Principal Component Analysis of all crocin data allowed to discriminate samples based on origin (Italy vs. other countries) and age and to confirm the feasibility of the use of crocin pattern as marker of quality of saffron.

Conclusions

A simple and fast procedure for the simultaneous determination of crocins has been developed by UHPLC-MS/MS. Multivariate analysis was used to distinguish saffron samples with different storage and process and data also confirmed that geographical traceability can be obtained by crocins concentration and isomers pattern.

Novel Aspect

The experimental data confirmed the feasibility to use -cis and -trans isomers of crocins as markers of quality, process and traceability.

References

1. Maggi L., Carmona M., Zalacain A., Kanakis C. D., Anastasaki E., Tarantilis P. A., Polissiou M. G., Alonso G. L., Food Research International, 43:1329-1334 (2010).
2. Melnyk J. P., Wang S., Marcone M. F., Food Research International, 43:1981-1989 (2010).

WP-87 / REAL-TIME AUTHENTICATION OF FOOD AND BEVERAGES USING DART-QDA LIVEID ANALYSIS

Sara Stead⁽¹⁾ - **Jayne Kirk**⁽²⁾ - **Kari Organtini**⁽³⁾

Waters Corporation, Scientific Operations, Manchester (1) - Waters Corporation, Product Management, Manchester (2) - Waters corporation, Scientific Operations, Milford (3)

Keywords: DART; mass detection; food fraud; chemometrics

Introduction

(Real-Time Authentication of food and beverages Using DART-QDa LiveID Analysis)

Currently there's no globally harmonised definition for "food fraud". Food and beverages of high commercial value are subject to fraudulent practice. Whisky is sold as the product of one distillery or as a blend. Analytical methods are required for process control and formulation. Whiskey characteristics are influenced by the cereals used in fermentation, distillation, maturation and blending regimes. Palm oil is also a commodity of socioeconomic importance.

Methods

DART was coupled to a single quadrupole MS to produce profile data. Linear X rail module was used for direct sample introduction. The DART source was operated with a helium gas flow of 1.5 L/min at 400°C. The MS was operated in full scan mode (mass range m/z 150-950) in positive and negative polarity; cone voltage 5-30 V; scan rate of 2 Hz. Chemometric models were generated using authentic samples, these models were used for real-time classification of whisky and palm oils.

Results

Within this study we have investigated the potential for the implementation of DART QDa LiveID as solution for the rapid authentication of high value food and beverage samples of both liquid and solid types with minimal sample preparation. Multivariate models for product authenticity have been generated using a collection of authentic reference samples of single malt whiskey and palm oils from different geographical locations. In whisky, ions with m/z between 50–300 we seen to be responsible for the brand level differences whereas ions at between m/z 300-400 drove the differences between the authentic and adulterated samples. In the palm oil, ionized species suspected to be generated from saturated fatty acids; carotenoids; vitamin A; tocopherols; tocotrienols; diglycerides and triglycerides were evident in the full scan spectral data. On the basis of these spectral features, it is possible to generate a multivariate model using a combination of PCA for data dimension reduction and LDA for class discrimination to distinguish between grades of palm oil processing and production type.

Conclusions

DART produces relatively simple mass spectra characterized by M^+ , $[M+H]^+$ in positive mode and M^- , $[M+H]^-$ in negative mode. Fragment ions are observed for some compounds, the degree of fragmentation can be influenced by the choice of gas. In both the whiskey and palm oil applications, positive polarity ionization was found to give most abundant spectral features and considered to be the more diagnostic mode for model training.

Novel Aspect

DART QDa with LiveID for rapid characterization and authentication of solid and liquid food and beverages

WP-88 / EVALUATION OF SOLID-PHASE MICROEXTRACTION FOR THE ISOTOPIC ANALYSIS OF VOLATILE COMPOUNDS IN APPLE

Lidija Strojnik⁽¹⁾ - *Federica Camin*⁽²⁾ - *Luca Ziller*⁽²⁾ - *Nives Ogrinc*⁽¹⁾

Jožef Stefan Institute, Department of Environmental Sciences, Ljubljana (1) - Fondazione Edmund Mach, Department of Food Quality and Nutrition, Research and Innovation Centre, San Michele all'Adige (2)

Keywords: SPME; GC-C-IRMS; apple; flavor compounds; authenticity

Introduction

The use of solid-phase microextraction (SPME) coupled with isotope ratio mass spectrometry (IRMS) for the $\delta^{13}\text{C}$ and $\delta^2\text{H}$ analysis of flavor compounds from apple fruit and recovery aromas has been evaluated using headspace sampling (HS). The method can be used to differentiate between natural and synthetic apple aroma compounds [1] and it can be easily transferred also to other commodities.

Methods

Factors affecting SPME were optimized [2], and the suitability of coupling SPME and gas chromatography-combustion-IRMS (GC-C-IRMS) [3] for determining $\delta^{13}\text{C}$ and $\delta^2\text{H}$ values was investigated. To do this, we investigated if HS-SPME extraction and analysis conditions result in isotopic fractionation. This was done by making a comparison of $\delta^{13}\text{C}$ and $\delta^2\text{H}$ values obtained by HS-SPME by GC-C-IRMS and $\delta^{13}\text{C}$, and $\delta^2\text{H}$ values determined using bulk analysis of reference material.

Results

Volatile components were extracted in 10 mL SPME vials containing 1 mL of the sample prior to gas chromatographic analyses. Optimal parameters were selected based on DVB/CAR/PDMS coated fibre and headspace sampling mode. Optimized SPME conditions were based on an adsorption time of 20 min at 30°C and desorption time of 60 s at 250°C. It is also shown that neither the analyte concentration nor the split/splitless mode affects $\delta^{13}\text{C}$ values. Finally, we verified that there are no matrix effects from the apple fruit or apple juice medium and that the technological production from apples to apple recovery aromas cause no isotopic fractionation. Therefore, comparison of $\delta^{13}\text{C}$ values is possible, and the use of cheaper apple samples instead of expensive recovery aromas for creating an apple aroma compound $\delta^{13}\text{C}$ database is appropriate. Such a database is required to provide reliable results for determining the authenticity of flavorings.

Conclusions

Our results show that by using an optimized SPME method we can obtain highly reproducible results and that SPME is, therefore, an appropriate tool for the isotopic analysis of volatile apple compounds extracted from fruit, juice or recovery aroma. Since many different compounds with different concentration are present in one sample the selection of reference material and appropriate processing and interpretation of the results obtained is crucial.

Novel Aspect

We report for the first time how different HS-SPME extraction parameters and analysis conditions affect the $\delta^{13}\text{C}$ and $\delta^2\text{H}$ values of apple aroma compounds obtained by GC-C-IRMS.

References

1. Van Leeuwen, K. A., Prenzler, P. D., Ryan, D., Camin, F. *Comprehensive Reviews in Food Science and Food Safety*, 13, 814-837 (2014).
2. Vas, G., Vekey, K. *Journal of Mass Spectrometry*, 39, 233-254 (2004).
3. Schipilliti, L., Bonaccorsi, I. L., Mondello, L. *Flavour Fragr. J.*, 32, 85-91(2017).

WP-89 / AROMATIC COMPOUNDS PROFILE OF DOC PIEDMONTESE WINES BY GC-MS WITH HS-SPME

Elisabetta Bonometti^(1,2) - ***Roberto Rabezzana***⁽¹⁾ - ***Giuseppina Cerrato***^(1,3) - ***Lorenza Operti***⁽¹⁾ - ***Eliano Diana***⁽¹⁾ - ***Francesca Turco***⁽¹⁾ - ***Emanuele Priola***⁽¹⁾ - ***Alessia Giordana***⁽¹⁾

University of Turin, Department of Chemistry, Via P. Giuria 7, 10125 Torino, Italy (1) - University of Turin, CrisDi Interdipartimental Centre for Crystallography, Via P. Giuria 7, 10125 Torino, Italy (2) - University of Turin, NIS Center for nanostructured interfaces and surfaces, Via P. Giuria 7, 10125 Torino, Italy (3)

Keywords: HS-SPME, GC-MS, wine, volatile organic elements, aromatic molecules

Introduction

Gas chromatography-mass spectrometry is widely employed to obtain important information on the aromatic profiles of wines, which are widely employed for differentiation and classification of wines based on geographical origin and grape cultivar. Wine is recognized to be very important in the human diet: the worldwide consumption of red wines reached over 240 mhl in 2017 [1].

Methods

On the basis of the general interest for local and sustainable products, we decided to focus our attention on some not yet studied Piedmont hillock wines. The volatile molecules were extracted by means of headspace microextraction (HS-SPME) with a triphasic fiber (CAR/PDMS/DVB 30/50 µm), and were quantified by means of gas chromatography-mass spectrometry. This extraction technique has gained universal approval and is employed in many applications such as aromas and fragrances recognition in food area, toxicology, environmental and biological matrixes [2].

Results and Conclusions

Multivariate chemometric approach allowed the identification of the representative profile of the selected typical Piedmontese wines, which will be useful for traceability and safeguard from food frauds.

References

1. OIV – Organisation Internationale de la Vigne et du Vin, april 2017 report
www.oiv.int/public/medias/5287/oiv-noteconjmars2017-en.pdf
2. S. Camara, M. Arminda Alves, J. C. Marques, Anal. Chim. Acta, 555, 191-200 (2006).

WP-90 / PROTEOMICS APPROACH FOR THE DIFFERENTIATION OF TUNA SPECIES

Ina Brümmer - ***Jens Brockmeyer***

University of Stuttgart, Food Chemistry, Stuttgart, Germany

Keywords: Authenticity, Food, Proteomics, HPLC-MS/MS, Tuna

Introduction

Worldwide studies prove that food fraud of fishes is a relevant issue[1][2]. Tuna species (*Thunnus* spp.) are often affected caused by limited availability and high price differences [3]. As methods for the differentiation of tuna species like polymerase chain reaction (PCR) and isoelectric focusing (IEF) are not sufficient, the aim of this project is to develop a method for tuna differentiation based on species-specific marker peptides.

Methods

Tryptical digested peptides of each species were measured with HPLC-MS/MS by information dependent acquisition (IDA). The identification of de novo peptides was done using the software PEAKS, comparison

was done manually. The identified potential marker peptides were verified by targeted experiments using HPLC-MS/MS. Finally, a targeted method for the differentiation of all relevant tuna species was developed.

Results

In this project, we included all eight tuna species (*Thunnus* spp.) and two close related tuna species, skipjack tuna (*Katsuwonus pelamis*) and mackerel tuna (*Euthynnus affinis*). After manually comparison of all de novo peptides, species-specific marker peptides which are present in all analyzed fishes of these species were found for all relevant species. After verifying these potential marker peptides, it was possible to develop a targeted method using HPLC-MS/MS.

In addition, marker peptides for the group of all three Bluefin tunas (*Thunnus thynnus*, *Thunnus maccoyii* and *Thunnus orientalis*) and also for the genus *Thunnus* were found.

Conclusions

Although tuna species are close related it had been shown that tuna species differentiation is possible using marker peptides analyzed with a targeted approach. In contrast to molecular biological methods, species mixtures can be also detected. After validation, this method could be implemented for routine analysis.

Novel Aspect

This developed method is the first method for the differentiation of all tuna species based on mass spectrometry using a Proteomics approach.

References

1. Espiñeira M. et al., *Journal of Agricultural and Food Chemistry*, 56, 10594-10599 (2008).
2. Jacquet J. L., Pauly D., *Marine Policy* 32, 309-318 (2008).
3. Kappel K., Schröder U., *Journal of Consumer Protection and Food Safety*, 10, suppl. 1., 31-34 (2015).

WP-91 / IDENTIFICATION AND QUANTIFICATION OF POLYPHENOL IN COCOA EXTRACT, STUDY OF THE EFFECTS OF FERMENTATION.

Guillaume Cazals - François Guéhi Gbogbri - S. Munier - Renaud Boulanger - Alain Michel - Christine Enjalbal - Patrick Poucheret

LMP, Université de Montpellier, Montpellier, France

Keywords: polyphenol, identification, quantification, cocoa, fermentation

Introduction

Regular consumption of chocolate is associated with beneficial effects on health [1]. Cocoa processing processes degrade the polyphenols it contains; the fermentation alone reduces by 90% the content of the polyphenols. It is therefore important to know its impact on the content and / or composition of polyphenols. Polyphenols were extracted from the post-harvest cocoa beans and the chemical compounds present in the cocoa extracts were identified and quantified by LC-MS

Methods

Identification of polyphenol were done by using UPLC separation coupled with high-resolution MSMS (QToF analyser). Identification of compounds were based on their exact mass and MSMS profile based on literature. Some of the flavanols identified were quantified by UPLC-MRM (triple quadrupole) with external calibration [2].

Results

The post-harvest cocoa beans were divided into two batches (fermented and unfermented). In order to study difference between these two batches we achieved to build a robust chromatographic method [3] in

order to get a good separation of compounds contained in the extracts.

In unfermented extract we observed twenty different compounds (e.g. Catechin, procyanidin B, ...) while in fermented extract we only observe seven compounds.

From these observed compounds we choose three of them (catechin, epicatechin, procyanidin B) well identified to achieve quantification in both extract. Quantification results show decrease in the concentration of this compound in fermented extract. For example, epicatechin is present at 5.9 µg/g of cocoa in unfermented extract and at 0.21 µg/g of cocoa in fermented extract.

Conclusions

Cocoa is a natural substance rich in polyphenols [4,5,6]. Fermentation is one of the essential processes in the processing of cocoa beans. This study demonstrated that this process degrades the polyphenols responsible for the health effects of cocoa. However, because of conflicting data reported in the literature, in-depth studies to fully identify different methods of cocoa fermentation are required [7]. In addition, to gain benefit, the research should take into account all the possible parameters such as the temperature, the fermentation time and the room in which the beans that deposit the initial content of the polyphenols during the fermentation processes are deposited. This would make it possible to avoid uncontrolled phenomena of natural, artificial or accidental origin that interfere during the fermentation of cocoa.

Novel Aspect

Relationship between post-harvest treatment of cocoa and its composition in polyphenol

References

1. Abbe M. M. J., Amin I., *Molecules*, 13, 2190-2219 (2008).
2. Afoakwa, E. O., Quao, J., Takrama, F. S., Budu, A. S.m Saalia, F. K., *International Food Research Journal*, 19, 1070-1077 (2012).
3. Andrés-Lacueva, C. ; Lamuela-Raventós, R. M. ; Jáuregui, O. ; Casals, I. ; Izquierdo-Pulido, M. ; Permanyer, J., *LC GC Eur*, 12, 902-904 (2000).
4. Albertini B., Schoubben A., Guarnaccia D., Pinelli F., Della Vecchia M., Ricci M., Di Renzo G. C., Blasi P., *Journal of agricultural and food chemistry* 63, 9948-9953 (2015).
5. Camu N., De Winter T., Addo S. K., Takrama J. S., Bernaert H., de Vuyst L., *J. Sci. Food Agric.* 88, 2288–2297 (2008).
6. Cienfuegos-Jovellanos E., Quinones M. D., Mu Steinberg erza B., Moulay L., Miguel M. Aleixandre A, 57, 6156–6162 (2009).
7. Counet C., Ouwerx C., Rosoux D., Collin S., *Journal of agricultural and food chemistry*, 52, 6243–6249 (2004)

WP-92 / INVESTIGATION INTO THE CIS ISOMERS OF DICHLOROGENIC ACIDS IN COFFEE

Silvia Colombari - Elena Guercia - Elisabetta De Angelis - Luciano Navarini

illycaffè spa, Aromalab, Trieste, Italy

Keywords: chlorogenic acids, isomers, UHPLC-MS, green coffee

Introduction

Chlorogenic acids (CGAs) are a large class of esters formed between quinic acid and hydroxycinnamic acids frequently present in plants as a complex mixture of positional and geometric (cis/trans) isomers, however not easily and univocally identified [1]. Additionally, the UV light exposure is proven to convert the trans double bond to cis and structure and activity relationship for these cis/trans isomers have not been the subject of detailed investigations in spite of the need to deepen the knowledge on these biologically active polyphenols [2,3].

Methods

Separation of analytes was achieved by UHPLC system equipped with a C18 column coupled to a Triple Quad MS, and, in order to discriminate isomers for dicaffeoyl quinic acids (diCQA), a targeted MS3 mode (m/z 515.0 + 353.0 on negative mode) experiment was used.

Ad hoc UV light exposure experiments were performed on standard solutions and green coffee extracts.

Results

Results confirmed that geometrical isomerism does not influence the fragmentation pattern, as previously reported in the literature [2], and an identification of the positional isomers is possible by means of an accurate chromatographic separation. Different behavior of light exposed green coffee extracts were evaluated, showing how differences in chlorogenic acids profile of different species could lead to different mixture of isomers after UV treatment.

Conclusions

This work showed that UV irradiation of diCQA samples produce three positional isomers, and separation of this isomers could be efficient with the use of a LC-MS/MS method. Differences in chlorogenic acids profile of various green coffee samples were showed and these findings could represent an important step in the detection and discrimination of isomeric plant metabolites.

Novel Aspect

UV irradiation experiment on green coffee extracts of different species with characterization of isomers of di-acylchlorogenic acids.

References

1. Clifford MN, Chlorogenic acids. In: Coffee. Vol 1. Chemistry. Clarke RJ; Macrae R (eds) Elsevier Applied Science Publications, London, UK (1985).
2. Clifford, N., Kirkpatrick, Jo., Kuhnert, N., Roozendaal, H., Food Chemistry 106: 379-385 (2008).
3. Baker E. S. et al, J Phys Chem Lett. 8, 1381–1388 (2017).

WP-93 / COMPARATIVE PROTEOMIC ANALYSIS OF THE METABOLIC FRACTIONS COMPOSITION IN MODERN AND OLD WHEAT GENOTYPES

Antonella Di Francesco⁽¹⁾ - ***Rosaria Saletti***⁽¹⁾ - ***Vincenzo Cunsolo***⁽¹⁾ - ***Vera Muccilli***⁽¹⁾ - ***Pasquale De Vita***⁽²⁾ - ***Salvatore Foti***⁽¹⁾

University of Catania, Department of Chemical Sciences, Catania, Italy (1) - CREA Cereal Research Centre, CREA Cereal Research Centre, Foggia, Italy (2)

Keywords: wheat, old genotypes, label free quantification, high resolution mass spectrometry

Introduction

Most of the modern wheat genotypes are derived from old wheats and have appreciable properties in terms of grain yield [1]. Wheat derived products plays an important role in many adverse reactions, such as celiac disease (CD), non-celiac wheat sensitivity (NCWS) and allergies. This has led to an increasing interest for the old wheat genotypes which are generally considered better tolerated than the modern ones, but without any scientific evidence.

Methods

Metabolic fractions were extracted from flours using 0.4 M NaCl, 0.067 M K₂HPO₄, pH 7.6 for 15 at 20°C and a final centrifugation at 12,000 x g for 15 min. The extraction was repeated twice. 50 µg of each protein extract were reduced, alkylated, digested and analysed by high resolution mass spectrometry. MS data

were processed by PEAKSQ and searched against a protein database, including only entries of Triticum, Oryza, Hordeum, Avena, Secale, Maize and Brachypodium.

Results

The aim of the present work is the qualitative and quantitative comparison by a proteomic approach of the metabolic protein fractions extracted from the mature kernel of two old Sicilian durum wheats landraces (Russello and Timiliarestre bianche), and Simeto, an improved durum wheat variety, widely spread in Italy and other Mediterranean countries, chosen as representative of the most widely commercial cultivars. The qualitative comparison of the protein composition revealed a remarkable similarity between old and modern cultivars. In particular, for each cultivar about 450 proteins were identified. The quantitative evaluation of the identified proteins detected that some protein including avenin, NLTPs and β amylase, are over-expressed in the old genotypes compared to Simeto. On the other hand, differential expression is also observed in the comparison of the two old cultivars. In particular α -amylase inhibitor 0.28 is over expressed in Timilia whereas serpin Z2A is over expressed in Russello.

Conclusions

This study represents the first exhaustive qualitative and quantitative characterization of the metabolic protein fractions of the old Sicilian landraces Russello and Timiliarestre bianche and its comparison with a modern cultivar.

Novel Aspect

This comparison would help to understand the relationship occurring between protein profile of old wheat varieties and potential benefits for human consumption.

References

1. Shewry P. R., J. Exp. Bot. 60, 1537-53 (2009)

WP-94 / PROTEOMICS GOING NUTS? PROTEOME PROFILING OF WALNUT KERNELS

Benjamin Dreyer - Dennis Kroesser - Hartmut Schlüter

University Medical Center Hamburg-Eppendorf, Mass spectrometric proteomics, Hamburg, Germany

Keywords: Food authenticity, Walnut, food fraud, reference protein barcodes

Introduction

Food authenticity is nowadays highly relevant since food fraud of expansive foods is a common problem. Differential proteomics approaches are powerful tools in this regard. Proteomics needs a tailored sample preparation for every kind of food. Especially foods with high content of lipids like nuts are challenging. In this study different sample preparation methods suitable for downstream LC-MS/MS analysis were tested for obtaining proteome profiles of the walnut.

Methods

In the 1st approach nutmeat was grinded and defatted with hexane. Proteins were extracted by boiling in a sodium deoxycholate (SDC) containing buffer and then proteolytically digested with trypsin. In the 2nd approach grinded nutmeat was boiled in SDC containing buffer for protein extraction. Proteins were precipitated with acetonitrile (ACN) and then digested with trypsin. Peptides were analyzed by LC-MS/MS and raw data was processed with Proteome Discoverer 2.0.

Results

The two protein extraction methodologies were compared regarding the quantity of the extracted protein

amount in relation to the weighed portion of nutmeat and the number of unambiguously identified proteins and peptides respectively. Protein amounts were determined with a BCA assay and the recovery per weighed in mg of nutmeat was calculated. The protein recovery was approx. 38 % higher if the nutmeat was defatted first. Nevertheless the protein yields of both methods were sufficient for LC-MS/MS analysis. Although the ACN precipitation procedure yielded less protein, it performed superior in terms of the number of peptide and protein identification. Approximately 3.5 times more peptides and about six times more proteins were identified by the 2nd approach using the ACN precipitation. Most abundant in all samples were seed storage proteins like legumin B-like and 11S globulin seed storage protein 2-like with up to 1000 peptide spectrum matches. The data suggests that the walnut nutmeat proteome has a large dynamic range.

Conclusions

Using ACN for precipitation of extracted proteins shows advantages over first extracting lipids with hexane for walnut sample preparation for proteomics applications. To overcome problems associated with large dynamic range proteomes 2-D LC approaches will be used to build spectral peptide libraries in the future and data independent acquisition will be used for reliable quantification to generate reference protein barcodes as fingerprints for authentic foods.

Novel Aspect

This study presents evidence that defatting of walnut nutmeat prior to protein extraction isn't beneficial for differential proteomics applications.

WP-95 / UNTARGETED HIGH RESOLUTION MASS SPECTROMETRY METHODS TO ASSESS SALMON AUTHENTICITY

Giuseppina Maria Fiorino⁽¹⁾ - *Marion Fresch*⁽²⁾ - *Ina Brümmer*⁽²⁾ - *Jens Brockmeyer*⁽²⁾ - *Ilario Losito*^(1,3) - *Linda Monaci*⁽¹⁾

National Research Council of Italy (CNR), Institute of Sciences of Food Production (ISPA), Bari, Italy (1) - University of Stuttgart, Institute of Biochemistry and Technical Biochemistry, Department of Food Chemistry, Stuttgart, Germany (2) - Department of Chemistry and SMART Inter-department Research Center, University of Bari "Aldo Moro", Bari, Italy (3)

Keywords: Salmon authenticity; Direct analysis in Real Time; Liquid Chromatography; High Resolution Mass Spectrometry; Multivariate statistical analysis.

Introduction

Two different untargeted approaches based on Direct Analysis in Real Time (DART) and Liquid Chromatography coupled to High Resolution Mass Spectrometry (LC-HRMS) were applied to the lipid and protein fraction of 50 farmed or wild-type salmons belonging to *Salmo salar* species to assess the feasibility of MS-based untargeted analysis for the discrimination between the two groups of salmons.

Methods

Salmon lipid extracts were analyzed by DART coupled to an Orbitrap mass spectrometer. For the proteomics approach, proteins were extracted using Urea buffer systems and subjected to trypsin digestion and sample cleanup. Samples of farmed and wild-type salmon were analyzed by LC-MS/MS on TT6600 triple TOF-MS in DIA and SWATH mode and data evaluation was performed employing different chemometric approaches.

Results

Lipids or peptides retrieved from DART-HRMS and LC-HRMS spectra, respectively, and exhibiting significant

spectral intensity were selected as variables for multivariate statistical analysis. The latter allowed to distinguish clearly wild-type from farmed salmons on the scatterplot based on the first two principal components, accounting for the 77% of the total variance.

The subsequent search in the online databases showed that some of the selected variables might belong to the fatty acids class. In particular, three saturated fatty acids (14:0, 16:0, 18:0) had the main discriminating loading for wild-type salmons, whereas unsaturated fatty acids (18:1, 18:2, 18:3) and several oxidized forms arising from them were found to have a higher incidence in farmed salmons, due to their relevant incidence in the feeding mixtures adopted in salmon aquaculture plants [1]. Similarly, samples were clearly distinguishable using proteomics approaches.

Conclusions

The application of untargeted DART/LC-HRMS approaches, integrated by multivariate statistical analysis, appears to be a promising tool for the discrimination of farmed and wild type salmons. DART/LC-HRMS spectra provide a set of variables able to distinguish salmons on the basis of their different living conditions.

Novel Aspect

A novel coupling between DART and Mass Spectrometry has been investigated and exploited in the untargeted mode to discriminate between wild-type and farmed salmons.

References

1. Blanchet C., Lucas M., Julien P., Morin R., Gingras S., Dewailly É., *Lipids* 40, 529 (2005).

Acknowledgements

The present research has received funding from the European Union's Seventh Framework Programme for research, technological development and demonstration, under grant agreement No. 613688 2 "Food Integrity".

WP-96 / CHARACTERIZATION AND QUANTITATION OF WATER-SOLUBLE VITAMINS IN WINE BY LIQUID CHROMATOGRAPHY-HIGH RESOLUTION AND TANDEM MASS SPECTROMETRY

Giovanni Inzalaco⁽¹⁾ - *Francesca Borghini*⁽²⁾ - *Alessandro Cavaglioni*⁽²⁾ - *Gianluca Giorgi*⁽¹⁾

Università degli Studi di Siena, Dipartimento Biotecnologie, Chimica e Farmacia, Siena, Italy (1) - ISVEA srl, Poggibonsi, Italy (2)

Keywords: water-soluble vitamins, wine, HPLC-HRMS

Introduction

Wine is the oldest alcoholic beverage obtained from the fermentation of grapes. Among micronutrients, vitamins play a role in wine production. Some have origin from grapes and/or are formed during fermentation, others maybe added for fine-tuning some properties or for avoiding undesired processes.

In this study we set up a simple and reliable method by using liquid chromatography-high resolution (HR) and tandem mass spectrometry (MS) for characterization and quantification of water-soluble vitamins (vitamins B, C and J) in Italian red and white wines. Furthermore, the method is applicable to beers, grapes and grape juices.

Methods

Preliminary studies by direct infusion have been carried out on vitamin standards using both HRMS (Q-Exactive Hybrid quadrupole-Orbitrap, Thermo Scientific) and an ion trap instrument (LCQ Deca, Thermo Scientific) both equipped with an ESI source in positive mode. MS/MS experiments by using higher-energy

collisional dissociation (HCD) in the Orbitrap and low-energy collision induced dissociations (CID) in the ion trap have been carried out to choose the appropriate reactions to be monitored. Successively, real samples have been filtered and injected in the HPLC equipped with a C18 column (Acquity BEH 100 × 2.1mm 1.7µm) and analyzed by HRMS. Analysis have been performed in Full MS-ddMS2; quantification by external standard method.

Results

Vitamin fragmentation patterns show differences between HCD and CID spectra. The former, acquired in high resolution, generally show a higher number of product ions and ions that in the CID spectra are shown only in MS3. The quantities measured in wine are not relevant in terms of daily vitamin intake, but the study is of oenological interest. Vitamin B12 and piridoxal were not detected and folic acid was found in trace only in beer samples. Ascorbic acid was detected only in white wines. Vitamin C is a quick-acting and powerful antioxidant that can be added in white winemaking to protect against light and short aerations, preserving freshness and fruity [1]. Also, thiamine is an allowed additive. Some compounds (*e.g.* pantotenic acid) were in higher amounts in grapes than in wines, indicating that these vitamins are used by yeasts [2]. Others (vitamin J, B2) showed a reverse trend probably because are not used by yeast alcoholic fermentation [3]and they are released in the medium as the result of yeast cellular lysis. Some are present only in wines (nicotidamide) and are produced by yeast metabolism during fermentation [4,5].

Conclusions

We developed an HPLC-HR-MS/MS method for the determination of water-soluble vitamins in wine and also can be extended to beer, grapes and grape juices. Vitamin concentrations are comparable with those reported in wines by other Authors. Different vitamin profiles were found in the analyzed matrices. Changes between grapes and wines might be attributed to the fermentation process.

Novel Aspect

Quite large-scale screening of water-soluble vitamins in Italian wines. Method fast, sensitive, specific and robust without any sample pretreatment.

References

1. Skouroumounis G. K., Kwiatkowski M. J., Francis I. L., Oakey H., Capone D. L., Peng Z., Duncan B., Sefton M. A., Waters E., J. Aust. J. Grape Wine Res. 11, 355-368 (2005).
2. Duc C., Pradal M., Sanchez I., Noble J., Tesnière C., Blondin B., PLoS ONE 12, 1-22 (2017).
3. Hucker B., Wakeling L., Vriesekoop F., J. Inst. Brew. 122, 126–137 (2016).
4. Orozco H., Matallana E., Aranda A., Mech. Ageing Dev. 133, 348–358 (2012).
5. Bitterman, K. J., Anderson R. M., Cohen H. Y., Latorre-Esteves M., Sinclair D. A., J. Biol. Chem. 277, 45099–45107 (2002).

WP-97 / SAMPLE PROFILER – A CLOUD-BASED APPLICATION FOR CONFIRMING FOOD AUTHENTICITY

Srividya Kailasam⁽¹⁾ - **Ralf Tautenhahn**⁽²⁾ - **Tim Stratton**⁽³⁾

ThermoFisher Scientific, LSMS, Bangalore, India (1) - ThermoFisher Scientific, LSMS, San Jose, United States

(2) - ThermoFisher Scientific, LSMS, San Jose, United States (3)

Keywords: Food-authenticity, Cloud-based application, Sample Profiler

Introduction

Confirming food authenticity is important not only for consumers and the food industry but also for government regulators. There is often the need to detect trace levels of contaminants in complex matrices, which requires the use of state-of-the art mass spectrometers. Sophisticated software tools are needed to

quickly process large data sets generated in these studies.

Thermo Scientific™ Sample Profiler is a new cloud-based software application that can be applied in diverse fields such as food authenticity, environmental monitoring, metabolomics, and detection of new designer drugs. Using this application, large data sets can be processed to detect and integrate compound peaks. Further, compound identification is automated by comparing the fragmentation spectra of compounds against mzCloud™, a web-based spectral database.

Using the Sample Profiler application, chemical profiles of samples of known quality and authenticity can be created within minutes. Individual samples of suspect/questionable quality can be readily compared against these profiles. Statistical evaluation of these comparisons helps users identify certain compounds of interest or 'marker compounds' that are present in significantly different concentration levels in the two sets – authentic and suspect samples. Thus, the presence or absence of marker compounds helps with the quick verification of authenticity and/or characterization of challenge samples when compared against profiles of samples of known/desired characteristics.

In this study, we used Sample Profiler to distinguish between soy sauce samples that were brewed during manufacturing from those that were not.

Methods

Sample Preparation

Different brands of soy sauces were purchased locally. Samples were prepared in triplicate by diluting aliquots of soy sauce with water (3:1). Pooled quality control (QC) samples were prepared by combining a fixed amount of each soy sauce into a pool and preparing a set of dilution samples as normal.

Test Method

Analysis was performed by LC-MS with separation achieved by reversed-phase linear gradient (ACN/ water, 0.1% formic acid) with a run time of 10 min. Mass spectrometric analysis was performed on a Thermo Scientific™ Q Exactive™ HF hybrid quadrupole-Orbitrap mass spectrometer operated in positive mode with electrospray ionization. Full scan MS1 (120,000 resolution FWHM at m/z 200) and precursor selected MS/MS (30,000 resolution) data was acquired. Sample injection order was randomized with QC samples injected every 10 injections.

Data Analysis

Mass spectral data for different types of soy sauces, such as those manufactured in different regions, stored in different types of containers, and even processed differently (by brewing or by not brewing), were collected. Data for blanks and QCs were also obtained. Raw data files were uploaded to Thermo Fisher Cloud and imported into the Sample Profiler application for further processing. The components were detected in samples after the application aligned all the retention times in all samples. Compound identification was done by matching fragmentation (MS2) spectra against the mzCloud spectral database. Sample Profiler aggregated the name, predicted formulae, m/z , RT, and area counts of each compound from all the samples.

Results

Background

Soy sauce is a popular condiment manufactured and consumed in many Asian countries. Traditionally, soy sauces are prepared by two methods: a) by naturally brewing fermented soya beans with wheat and brine, which can take several months, or b) by using hydrolyzed vegetable protein (HVP), which takes only days. Soy sauces manufactured by the two processes differ in their physical characteristics such as aroma, flavor, and texture. In addition, soy sauces made by acid-hydrolysis of soy proteins have a longer shelf life.

In this study, different brands of soy sauces purchased locally were analyzed to determine the differences in their compositions and identify spurious samples based on our knowledge of genuine products. HRAM data of each sample was collected in triplicate followed by statistical analysis using Sample Profiler, an application hosted in Thermo Fisher Cloud.

Data Processing in the Cloud

Large data sets uploaded to the Data Manager in the Thermo Fisher Cloud can be easily added to projects created in Sample Profiler and processed using suitable processing methods. Data processing involves

aligning the chromatographic peaks, followed by component extraction. The accurate mass information was used to predict elemental composition for detected components. Finally, compounds were identified by matching the predicted formula and MS2 fragmentation data against ChemSpider and mzCloud™ databases, respectively.

Creating Profiles

The application also allows the user to aggregate samples with similar characteristics such as packaging and or manufacturing process which can be pre-defined as “Study Factors” prior to data analysis. These sample subsets that share one or more study factors can be then used to establish a ‘product-profile’. Other samples having same or even different study factors can be compared against chosen profiles to evaluate the extent of similarity.

After all samples including QCs and blanks are processed, two databases of compounds or profiles, a) All Samples and b) All QCs samples, are automatically created by the application. The user can also create additional profiles or subsets of the All Samples profile by partitioning it based on study factors or criteria such as manufacturing region or packaging material used to characterize the samples. For example, in the present study, a database or a profile of compounds present in soy samples that were brewed during manufacture was created.

In the present study, a profile of soy sauces that were brewed during manufacturing was used to determine the quality of soy samples that were not brewed. The data indicated that the chemical compositions of both types of soy sauce are similar in that they contain same compounds in similar concentration range. Yet, some compounds were either not found at all or found in much lower concentrations in soy sauces that were not brewed. Similarly soy samples manufactured in different regions or samples stored in different containers can be compared and statistically distinguished from one another.

The application displays the profile of “Brewed” soy-sauces, in a tabular format which includes details such as compound name, predicted formula, MW, RT, and highest intensity value from all samples for 11,572 compounds found in this subset of soy sauces. The spread in the intensity values found across all samples is shown as distribution plots. The ‘Compounds table’ can be filtered to show the compound details in any of the profiles created in the project.

Trend line Plots

Sample Profiler also allows for visualization of the trends in compound areas. When trend lines for compound areas across all samples in the profile are plotted, those samples that may have significantly higher or lower concentrations of one or more compounds can be easily identified. These samples can then be investigated further to get a better understanding of the products’ chemical composition.

Comparing Different Types of Samples

Using Sample Profiler, comparisons between different types of soy sauces can be easily set up and data processed at different fold changes and adjusted p-values. In the present study, a profile of brewed soy sauces was created and used for comparing against not-brewed samples. Data analysis performed by setting the fold change value to 8 and adjusted p-value to > 0.001 . It was found that of a total of 11,803 compounds, the not-brewed samples were found to contain 267 compounds in significantly higher abundance, and 1791 compounds in significantly lower abundance than in brewed samples.

Marker Compounds

Once a product’s profile is well understood, certain compounds when present in a product can be labeled as markers of authenticity or contaminant. These marker compounds are not only useful in qualifying study samples but also for setting up targeted data acquisition.

In this study, Sample Profiler was used to search for genistein and glycitein specified by name. Post search, the table was updated with formulae, MW, and RTs of the two compounds along with the total number of samples and QCs in which they were found. After the compounds are added to the table, Sample Profiler will automatically search for their presence in any comparison set up in the project.

The relative distributions of the two marker compounds in the comparison between brewed and not-brewed soy sauces were determined. From the results obtained, brewed soy sauces are understood to have higher abundances of these two marker compounds.

This knowledge can be applied to quickly verify the authenticity of a challenge sample(s) by setting up appropriate comparisons between them and known profiles and studying the relative distributions of

marker compounds.

Conclusions

Sample Profiler helps in the creation of chemical profiles of authentic samples and identification of 'marker compounds'. Marker compounds help with easy verification or characterization of suspect/challenge. The application also allows for visualization of trends in the data in addition to TIC and XIC traces. In the present study, it has been shown that soy sauces manufactured in different regions or packaged differently can be easily differentiated by the Sample Profiler application. Compounds such as glycitein and genistein were found to be present in higher abundance in brewed soy sauces as compared to not-brewed samples.

Novel Aspect

Sample Profiler, a cloud-based data analysis software that can be process large data sets has been used for establishing food authenticity.

References

1. Ralf Tautenhahn, Rahel Comte , Tim Stratton, Srividya Kailasam, Dipankar Ghosh and Heinz Singer, Cloud-based analysis of complex sample systems for outlier and trend analysis in the context of water monitoring, ASMS Poster (2016)

WP-98 / RAPID AND DIRECT ANALYSIS OF COFFEE BEANS (GREEN AND ROASTED) AND ONLINE MONITORING OF COFFEE ROASTING BY SIFT-MS: FOCUS ON VOLATILE ALDEHYDES

Juliia Khomenko⁽¹⁾ - *Andrea Caretta*⁽²⁾ - *Valentina Lonzarich*⁽³⁾ - *Franco Biasioli*⁽¹⁾ - *Luciano Navarini*⁽³⁾

Edmund Mach Foundation, Department of Food Quality and Nutrition, San Michele all'Adige, Italy (1) - SRA Instruments S. p. A., Cernusco sul Naviglio, Italy (2) - Aromalab illycaffè S. p. A., Research & Innovation, Trieste, Italy (3)

Keywords: direct injection, volatile organic compounds

Introduction

Aroma of coffee beans (*Coffea arabica*) contributes to its quality evaluation and its origin tracing. The part of coffee aroma such as volatile aldehydes is linked directly to volatile organic compounds (VOCs) naturally occurring in green coffee beans as well as to roasting process. Fast and direct volatile aldehydes detection in green beans, during coffee roasting, and in roasted beans by a selected ion flow tube mass spectrometer (SIFT-MS) may possibly be a way for characterization of coffee origins.

Methods

VOCs of coffee beans of 5 geographical origins and 3 batches each were measured by SIFT-MS (SYFT VOICE200 ultra, Syft Ltd, New Zealand). 5 green coffee beans were placed in a 22 mL glass vial, incubated at 37°C for 30 min, and measured for 1 min by an autosampler coupled to SIFT-MS. Single bean was placed in a vial and roasted in the oven at 211°C for 20 min. SIFT-MS monitored roasting through the oven ventilation. Then vials were closed and measured in the same way as green beans.

Results

Twenty one aldehydes, two alcohols, and two methoxypyrazines were monitored in green beans along with a unspecific monoterpenes peak. Some typical compounds such as furans, pyridine, and pyrazine were added to the list of monitored compounds for the online monitoring of coffee bean roasting and the screening of roasted coffee beans. During roasting two types of evolution curve of VOC emissions were observed: the peak of emission at the beginning (2-nonenal, benzaldehyde, and 2-pentylfuran) and at the end of the roasting (hexanal, 2-hexenal, and others).

The special attention was paid for the possibility of aldehydes to discriminate the coffee bean origin. In

green beans origin was characterized by the significant differences in the concentrations of five aldehydes such as benzaldehyde, hexanal, 2-methyl-2-propanal, 2-methylbutanal, and pentanal. The concentration of aldehydes in roasted beans was significantly different from green ones but no significant difference according to their origin was observed.

Conclusions

The technique was successfully applied for screening of green and roasted coffee beans and for online monitoring of the coffee bean roasting. The main differences in aldehyde content in coffee beans were found in green beans rather than in roasted ones.

Novel Aspect

For the first time the analysis of green and roasted coffee beans and online monitoring of VOC emission during coffee roasting was performed by SIFT-MS. Moreover for the first time SIFT-MS was coupled with a multipurpose GC automatic sampler.

WP-99 / DIFFERENTIAL PROTEOMIC PROFILING OF TRUFFLES

Dennis Kroesser - Benjamin Dreyer - Hartmut Schlüter

University Medical Center Hamburg-Eppendorf, Mass Spectrometric Proteomics, Hamburg, Germany

Keywords: Proteomics, Truffle, LC-MS/MS, Food Authenticity

Introduction

Food fraud is becoming an increasingly critical topic. For expensive foods like truffle, incorrect labeling or contamination of premium products with inferior products are common. Truffles are scarcely researched and thus detection of food fraud challenging. Here, protein extraction methods were tested and a LC-MS/MS based differential proteomics approach was developed for determining reference LC-MS/MS data for future validation of truffle quality.

Methods

Homogenates from authentic European and Asian truffle samples were obtained by picosecond infrared-laser (PIRL) ablation or by lyophilization and grinding of truffles. Protein extracts were boiled in presence of sodium deoxycholate and enzymatically degraded by trypsin. Peptides were analyzed by LC-MS/MS in data-dependent acquisition mode to generate spectral libraries and, for obtaining reference LC-MS/MS data, in data-independent acquisition mode.

Results

The quantitative comparison of homogenates from PIRL ablation and classical homogenization of truffle showed an at least twofold change in relative amount for 15 % of extracted proteins. There was no preference towards a certain method. In the differential proteomic analysis of truffles by measurements in DIA mode, reference LC-MS/MS data were generated for the samples available at the moment. By hierarchical clustering analysis and principal component analysis samples from European truffles could be well distinguished from Asian truffles.

Conclusions

The quantitative differential proteomics approach generated reference LC-MS/MS data by which truffles of different origin can be distinguished. Therefore this approach is a future tool for food fraud detection.

Novel Aspect

This is the first study presenting LC-MS/MS-based data as reference for validation of truffle quality.

WP-100 / A NOVEL ROBUST PROBE FOR THE RAPID, ON-LINE CHARACTERISATION OF FOOD SAMPLES USING RAPID EVAPORATIVE IONISATION MASS SPECTROMETRY

Andras Marton⁽¹⁾ - *Julia Balog*⁽¹⁾ - *Richárd Schöffler*⁽¹⁾ - *Steven Pringle*⁽²⁾ - *Zoltán Takáts*⁽³⁾

Waters Research Center, Budapest, Hungary (1) - Waters Corporation, Wilmslow, United Kingdom (2) - Imperial College, London, United Kingdom (3)

Keywords: Ambient ionization method, Rapid Evaporative Ionization Mass Spectrometry (REIMS), Food authenticity,

Introduction

There has been a significant increase of the incidence rate of food fraud cases related to products originating from uncontrolled origins. Rapid Evaporative Ionization Mass Spectrometry (REIMS) has been used for the analysis of various meat samples and has shown to be capable of the fast identification of diverse tissue types based on their lipid fingerprints. We aimed to test the feasibility of REIMS combined with our robust sampling probe for different food products.

Methods

Authentic mangalica pork and carp fish samples from 6 different suppliers based in Hungary were sourced. All samples were analysed using a custom built bipolar sampling probe. The aerosol generated by the thermal ablation of the samples was introduced into the REIMS source of a Xevo G2-XS ToF MS and spectra were acquired between 150-1500 *m/z*. By the application of multivariate algorithms, our goal was to identify the species or location of origin of each sample.

Results

The original monopolar sampling probe was modified forming a robust bipolar probe to allow rapid on-line sampling of products. The spectra acquired from the samples were first pre-processed (lockmass correction and background subtraction), then subjected to Principal Component Analysis with the results processed using Linear Discriminant Analysis. All classification models were tested on leave-one-sample out cross-validation, and an independent validation set. In the case of the pork study, the cross-validation resulted in >80 % accuracy, with 12 out of 15 samples correctly identified from the independent validation set. The most significant peaks in the spectra were identified using tandem MSMS measurements and featured mainly fatty acids, glycerophospholipids and triglycerides. In case of the carp study, the location of origin showed good separation in the first three principal components, and the cross-validation resulted in an 85% correct classification rate.

Conclusions

The results demonstrate that the lipidomic profiles can be taken in the few seconds timeframe and the profiles show good species- and origin- level specificity. The custom built food probe was found to be suitable for on-line sampling in an industrial environment. REIMS technology was found to be a promising tool for food safety applications providing a reliable and simple method for the rapid characterization of food products.

Novel Aspect

On site, on-line rapid sample identification with robust sampling probe based on REIMS lipid fingerprints.

WP-101 / CORRELATING CONSUMER SENSORY EXPERIENCE OF APPLE AROMA WITH VOC PROFILES ACQUIRED BY AUTOMATED TD-GC-MS

Gareth Roberts⁽¹⁾ - *David Barden*⁽¹⁾ - *Nick Bukowski*⁽²⁾ - *Laura McGregor*⁽³⁾ - *Massimo Santoro*⁽¹⁾ - *Natasha D. Spadafora*⁽⁴⁾

Markes International, Marketing, Llantrisant, United Kingdom (1) - SepSolve Analytical, Product Development, Peterborough, United Kingdom (2) - SepSolve Analytical, Marketing, Peterborough, United Kingdom (3) - Università della Calabria, Markes International, Food Sciences, Rende (CS), Italy (4)

Keywords: automated thermal desorption, apple cultivars, sensory analysis, VOCs

Introduction

The aroma of fruit is regarded as one of the most important quality parameters sensed by consumers, and is influenced by numerous factors including genotype and growth conditions. In this study we determined the VOC profiles of apple cultivars and correlated them to consumer sensory evaluation.

Methods

Apple headspace (cultivars 'Smitten' and 'Granny Smith') was sampled onto sorbent tubes using a micro-chamber system, with analysis by thermal desorption–gas chromatography–mass spectrometry. Twenty panellists were asked to rate the samples according to nine sensory characteristics, and two-dimensional hierarchical cluster analysis(HCA)was performed to relate these scores to the VOC profile.

Results

A total of 29 VOCs were identified across the samples, of which 20 were esters. 'Smitten' was found to have significantly higher hedonic scores on all sensory measurements than 'Granny Smith', and significant positive correlations were found between high 'Taste/Flavour', 'Fruitiness' & 'Aroma/Smell' scores and five of the esters identified in the headspace analysis. These compounds are: (with their odour as described in an olfactometry study of apple cultivars¹) 2-methylpropyl acetate (pear, apple), n-butyl acetate (sweet, fruity), 2-methylbutyl acetate (fruity, apple), n-pentyl acetate (fruity, banana), and n-hexyl acetate (fruity, pear). In that study, the first three of these compounds were identified as 'principal odorants'.

Conclusions

We conclude that esters could be used as objective markers for the consumer appreciation of fresh apples, enabling optimum conditions for processing and storage to be identified without recourse to expensive sensory panels in every case. Such tests could also be used as part of routine quality control by the producer and retailer.

Novel Aspect

Detection of trace-level VOCs was only possible because of the large micro-chamber sample volumes and the low thermal desorption split ratios.

References

1. E. Mehinagic, G. Royer, R. Symoneaux, F. R. Jourjon, C. Prost, Characterization of odor-active volatiles in apples: Influence of cultivars and maturity stage, *Journal of Agricultural and Food Chemistry* 54: 2678–2687 (2006).

WP-102 / FOOD PROTEOMICS: A NOVEL APPROACH TO FOOD AUTHENTICATION USING GLOBAL PROTEIN SIGNATURES AND SPECIES SPECIFIC MARKER PEPTIDES

Laura Viet - Pauline Krimm - Jens Brockmeyer

University of Stuttgart, Food Chemistry, Stuttgart, Germany

Keywords: proteomics, food authenticity, LC-MS/MS, global marker peptides

Introduction

Food fraud by undeclared blending of meat from different species is a highly relevant topic in food

industries as shown by the European horse meat scandal in 2013 [1]. For the detection of fraudulent blending with undeclared species to the earliest stage highly efficient methods for identification and quantification are necessary.

Methods

Peptide based high performance liquid chromatography-multiple reaction monitoring methods using species specific marker peptides allow the detection of less than 1 % horse meat or pork in a beef matrix [2]. Additionally, unspecific marker peptides - present in numerous species - have been discovered. Both specific and unspecific marker peptides have been discovered using information dependent acquisition and subsequent de novo sequencing of tryptic peptides.

Results

In this study, we investigated the ratios of species specific markers and global markers for different binary mixtures of beef and pork, beef and horse and pork and horse as well as the amounts of the global markers in different mammalian meat products by means of high resolution-multiple reaction monitoring. Considering the constant ratio of specific and global markers for each species, these ratios can be used to detect alterations of food products. While the amount of species specific marker peptides decreases in case of food fraud, the amount of global markers remains more or less the same. Thus, a decrease of the ratio can be used to prove the presence of species mixtures even in cases where the substituting species are not known (untargeted food authentication). Furthermore, a linear correlation between this change of ratios and the extend of the alteration offers an approach to quantify the addition of undeclared meat. The results show that for most cases a good linear correlation with correlation coefficients of $R^2 > 0,95$ can be achieved.

Conclusions

Our data shows that the targeted analysis of species specific and global markers and the evaluation of their ratio is a promising way of detecting food fraud and additionally relatively quantifying fraudulent blending of meat from undeclared species.

Novel Aspect

This food authentication method uses unspecific markers instead of species specific ones allowing the detection of food fraud even when the substituting species are not known.

References

1. O'Mahony P. J., QJM, 106, 595-597 (2013).
2. von Bargaen C., Dojahn J., Waidelich D., Humpf H. U., Brockmeyer J., J. Agric. Food Chem, 61, 11986-11994 (2013).

WPS-S05 - MS tools in doping & toxicology

Chairs: Mario Thevis, Gianluca Bartolucci

WP-103 / DEVELOPMENT AND VALIDATION OF MULTI-TARGET ANALYTICAL PROCEDURES FOR THE IDENTIFICATION OF POTENTIAL DOPING AGENTS IN 'OVER THE COUNTER' PRODUCTS

Fabio Comunità⁽¹⁾ - **Francesco Botrè**^(1,2) - **Martina Capozzi**⁽¹⁾ - **Xavier de la Torre**⁽¹⁾ - **Monica Mazzarino**⁽¹⁾

Laboratorio Antidoping Federazione Medico Sportiva Italiana, Rome, Italy (1) - Dipartimento di Medicina Sperimentale, "Sapienza" Università di Roma, Viale Regina Elena 324, 00161 Rome, Italy (2)

Keywords: Multi-target procedures, Doping analysis, Nutritional supplements, Designer drugs

Introduction

The regular use of nutritional supplements by athletes is a very common practice [1]. Contamination and/or adulteration of these products has been claimed to be the cause of numerous adverse analytical findings [2]. Different research groups confirmed the unlabeled presence of anabolics and stimulants in various products [3-4]. For the above reasons it becomes necessary to develop methods for the detection of the widest variety of substances included in the WADA list [5].

Methods

Both liquid and gas chromatography coupled to mass spectrometry were utilized. Different pre-treatment techniques including liquid-liquid extraction and solid phase extraction were evaluated. Validation was evaluated through specificity, sensitivity, repeatability, and recovery.

Results

Two different analytical protocols were developed: the former was specifically designed for the analysis of anabolic substances extracted from nutritional supplements by liquid-liquid extraction with n-pentane, and the latter for the remaining compounds included in the Prohibited List by solid phase extraction with weak cationic exchange columns. The procedures were validated and presented high specificity (no interferences at the retention time of all analytes and no carryover), sensitivity (all compounds could be detected at or below 20 ng/g), repeatability (variability of the relative retention times lower than 0.1% and of relative abundances of selected ion transitions lower than 15% for all the analytes considered in the study) and good recoveries (over 60% for all analytes). The efficacy of the multi-analyte procedures here developed was tested by analysing different products. Results revealed the presence of testosterone in a dietary supplement, while higenamine, a β 2-agonist, was found in an herbal massage cream.

Conclusions

Two different multi-target analytical methods (the former based on GC-MS and the latter on LC-MS) were designed, developed, validated and successfully applied. Both procedures allow to detect the presence of doping compounds in different products, confirming that both the newly proposed protocols match the requirements of routine analytical methods for the detection of trace amounts of prohibited substances in products available “over the counter”.

Novel Aspect

Two complementary newly developed methods allow the detection of banned substances in both nutritional supplements and in “over the counter pharmaceutical” formulations.

References

1. Maughan R., Depiesse F., Geyer H., Journal of Sports Science, 25(S1), 103 (2007).
2. Maughan R., Journal of Sports Sciences, 23, 883 (2005).
3. Parr M. K., Pokrywka A., Kwiatkowska D., Schänzer W., Biology of Sport, 28, 153 (2011).
4. Perrenoud L., Saugy M., Saudan C., Journal of Chromatography B, 877, 3767 (2009).
5. WADA. The World Anti-Doping Code International Standard. Prohibited List, January 2018. Available at: https://www.wada-ama.org/sites/default/files/prohibited_list_2018_en.pdf.

WP-104 / LARGE VOLUME INJECTION: IMPROVING THE DETECTION CAPACITY OF ENDOGENOUS ANABOLIC STEROIDS BY GC-C-IRMS

Loredana Iannella⁽¹⁾ - ***Cristiana Colamonicì***⁽¹⁾ - ***Davide Curcio***⁽¹⁾ - ***Francesco Botrè***^(1, 2) - ***Xavier de la Torre***⁽¹⁾

Laboratorio Antidoping, Federazione Medico Sportiva Italiana, Largo Giulio Onesti 1, 00197, Rome, Italy (1)
 - Dipartimento di Medicina Sperimentale, “Sapienza” Università di Roma, Viale Regina Elena 324, 00161 Rome, Italy (2)

Keywords: Anabolic androgenic steroids, Doping control, Isotopic Ratio Mass Spectrometry, Programmed Temperature Vaporizer Injector

Introduction

The confirmation of the synthetic origin of endogenous anabolic androgenic steroids is performed in all doping laboratories accredited by WADA through a GC-c-IRMS analysis, due to the different carbon isotope abundance, $^{13}\text{C}/^{12}\text{C}$ ($\delta\text{‰}$), between endogenous steroids and synthetic ones. The high sensitivity requirements of current method and the limited volume of urine disposable for the analysis, led us to improve the injection conditions using the PTV technique [1].

Methods

Blank urine samples were spiked with boldenone and 19-norandrosterone standards and cleaned up by two sequential HPLC steps. The purified samples were injected in a GC-c-IRMS instrument, replacing the splitless injection mode with the use of a programmed temperature vaporizing injection (PTV). The initial sample volumes of urine was fixed at 6 mL, whereas different final injection volumes (2 to 10 μL) were evaluated.

Results

At present, the anti-doping laboratories must activate the confirmation analysis for boldenone and 19-norandrosterone by GC-c-IRMS whenever their concentration is greater than the limits defined in reference Technical Documents (TD). In details the "TD2016IRMS" [2] and the "TD2017NA" [3] impose the GC-c-IRMS confirmation for samples with a concentration in the range 5-30 ng/mL for boldenone and 2-15 ng/mL for 19-norandrosterone. In the ordinary method, 21 mL of urine are used, leading to a final injected volume, in splitless mode, of 2-3 μL [4,5]. These conditions are necessary to obtain adequate signals within the linear range of the instrument (400 to 7000 mV) considering the low LOD of boldenone and 19-norandrosterone (4 and 2 ng/mL, respectively). We investigated the possibility to increase the injection volumes, reducing progressively the initial urine volumes. We noticed that the compounds response signals improved proportionally to the larger sample volumes injected, without any significant increase of the background noise.

Conclusions

The results proved that the possibility to adjust the injection volumes, choosing injection volumes larger than those traditionally used in the splitless mode, allows to reduce and standardize the initial urine sample volumes, making faster and less laborious the preliminary sample preparation steps. The reduced urine volume also ensures a reduced risk of column overloading and a decrease the matrix effect in the instrumental analysis.

Novel Aspect

The PTV technique could be extended to detect the exogenous origin of other compounds with a very low concentration.

References

1. Polet M., Van Gansbeke W., Deventer K., Van Eenoo P., *Biomedical Chromatography*, 27, 259–266 (2012).
2. WADA Technical Document – TD2016IRMS Detection of Synthetic Forms of Endogenous Anabolic Androgenic Steroids by GC/C /IRMS.
3. WADA Technical Document – TD2017NA Harmonization of analysis and reporting of 19-Norsteroids related to nandrolone.
4. De la Torre X., Curcio D., Colamonici C., Molaioni F., Botrè F., *Drug Testing and Analysis*, 5, 834–842 (2013).
5. De la Torre X., Colamonici C., Curcio D., Molaioni F., Pizzardi M., Botrè F., *Steroids*, 76, 471–477 (2011).

WP-105 / IS GC-MS ANALYSIS SUFFICIENT FOR IDENTIFICATION OF ILLICIT DRUGS FROM FESTIVE GATHERINGS?

Sophie Magréault⁽¹⁾ - *Christel Grondin*⁽²⁾ - *Catherine Lamoureux*⁽³⁾ - *Emmanuel Bourgoigne*⁽⁴⁾

AP-HP-hopital Robert Debré, Laboratoire de Pharmacologie, Paris (1) - AP-HP-hopital Lariboisière, Laboratoire de Toxicologie, Paris (2) - Service Commun des Laboratoires, Laboratoire de Paris, Paris (3) - AP-HP-hopital Lariboisière/Université ParisDescartes, Laboratoire de Toxicologie, Paris (4)

Keywords: GC-MS; clinical toxicology; illicit drug; Festive drugs

Introduction

Poisoning by recreational use of psychoactive substances is increasingly common in Europe and is responsible for death. In toxicology, identification needs different extractions procedures and derivatization before GC-MS analysis. Developing systematic strategies for screening drug compounds in tablets, powders, herbs or liquids remains a challenge. Through different French festive events, samples were collected with the main objective to discuss advantages and limitations of such systematic analytical procedures using conventional GC-MS.

Methods

For solids, analytes are dissolved in EtOH. For liquids, analytes are dissolved in water and MeOH, pH was measured and acidic and basic liquid-liquid extractions were performed. After evaporation, samples were reconstituted in ethanol before injection onto a GC-MS operated in scan mode. Derivatization (acetylation, silylation, propylchloroformate) are performed and each spectra compared to different commercial (MPW, SWGDrugs, NIST) or in-house libraries.

Results

Between July and October 2016, 81 recreational substances were collected during 10 interventions in festive gatherings. Using a systematic workflow, drugs were found in every sample, whatever the matrix (solid, liquid, resin, oil) except in blotters where LSD was expected. Drugs found were mainly "classic drugs" such as cocaine, amphetamines, or cannabis. Different derivatization is mandatory for complete identification of drugs enabling increasing sensitivity for low concentration diluents, cutting agents, adulterants and other psychoactive substances. Propylchloroformate derivatization allows covalently binding on Nitrogen heteroatom, whereas acetylation or silylation bind on oxygen heteroatom. As an example, in "Cocaine" samples, levamisole, phenacetin, caffeine, diltiazem, hydroxyzine were detected. Three new psychoactive substances (NPS) were also found: ethylphenidate, dibutylone (bk-DMBDB) and thiopropamine (tieno-amphetamine) highlighting the need for up-to-date and different mass spectra libraries.

Conclusions

The applied decision tree allowed finding all except one drugs in the samples, whatever the matrices. Drawbacks are the hard laboratory work, choosing the dilution identifying major analytes but also low concentrated cutting agents; different derivatization for exhaustive identification. Up-to-date different mass spectra libraries are also mandatory for NPS drug discovery. Without the latter, skilled technicians are required. GC-MS is not devoid of limitation. A designer drug (1p-LSD) present in blotting paper could not be identified and accurate mass analysis solved this issue. Ambient pressure ionization coupled to accurate mass is definitely the present/future for direct sample analysis. Nevertheless, conventional and systematic approaches efficacy should not be underestimated.

Novel Aspect

Application of a systematic decisional tree for identification of illicit drugs in solid/liquid/oil/plants samples by GC-MS in clinical toxicology

WP-106 / EVALUATION OF THE METABOLIC PROFILE OF 2-METHIOPROPAMINE BY LC-MS/MS IN MICE AND IN HUMAN LIVER MICROSOMES

Cristian Camuto^(1,2) - ***Sheila Pellegrini***⁽²⁾ - ***Monica Mazzarino***⁽¹⁾ - ***Xavier de la Torre***⁽¹⁾ - ***Francesco Botrè***^(1,3)
- ***Matteo Marti***⁽⁴⁾ - ***Fabio De Giorgio***⁽⁵⁾

FMSI Antidoping, Rome, Italy (1) - University La Sapienza-Department of Chemistry, Rome, Italy (2) - University La Sapienza, Department of Experimental Medicine, Rome, Italy (3) - University of Ferrara, University of Ferrara-Department of Experimental Medicine and Surgery, Ferrara, Italy (4) - University La Cattolica, Institute of Public Health, Rome, Italy (5)

Keywords: 2-Methiopropamine, Drug Metabolism, Liquid Chromatography-Mass spectrometry, Mice

Introduction

The aim of the present study was to investigate the phase I and phase II metabolism of 2-methiopropamine, a structural analogue of methamphetamine [1-3], after acute administration in mice, in order to select the most appropriate marker(s) of intake, also defining the excretion windows in urine. In vitro metabolism studies were carried out in order to enzymatically synthesize the metabolites of MPA and to set up the most appropriate sample pre-treatment procedures for LC-QqQ analysis.

Methods

A dose of 10 mg/kg was selected for the in vivo metabolism studies of MPA and urine samples were collected every 3 hours in the range of 0-9 hours after the injection. The In vitro studies were carried out with HLM (human liver microsomes). Samples from both studies were analysed, after extraction with tert-butyl-methyl ether at pH 9, using an Agilent 1200 HPLC equipped with a SUPELCO C18 column coupled to an API4000 triple quadrupole (Sciex) with ESI source, operating in positive-ion mode.

Results

To reproduce the in vivo metabolism, the in vitro metabolism protocol was optimized considering different conditions: substrate concentration 20 µM, proteins concentration 0.5 mg/mL, phosphate buffer 0.1 M at pH 7.4, and an incubation time of 4 h at 37°C. Three main metabolites were identified: Nor-MPA, Hydroxy-MPA, and Nor-hydroxy-MPA. To optimize sample pre-treatment different extraction solvents (tert-butylmethyl ether, diethyl ether and ethylacetate) and pH values (7, 9 and 12) were evaluated. The best recoveries (higher than 70%) for the principal metabolites were obtained at pH 9. In vivo metabolism studies shown the formation of three principal phase I metabolites identified as: Nor-MPA, Hydroxy-MPA, and Oxo-MPA. 2-MPA and its metabolites show a maximum of excretion in the first 3 h from administration to mice showing an increasing conversion into phase II metabolites both glucurono and sulfo-conjugates for MPA, Nor-MPA, Hydroxy-MPA; Nor-Hydroxy-MPA was detected only as a phase II metabolite.

Conclusions

The excretion profile of 2-MPA after acute administration was evaluated identifying its principal metabolites and their excretion windows in mice. MPA and Nor-MPA are the best markers of intake in the first 0-24 hours after the administration of MPA. The present method was qualitatively validated for MPA and Nor-MPA according to WADA technical document[4-5] defining: recovery, LOD, selectivity and repeatability.

Novel Aspect

There are few studies about MPA metabolism [3]. This is the first study identifying MPA excretion profile including phase I and phase II, also selecting the best markers of intake.

References

1. Blicke FF, Burckhalter JH. Journal American Chemistry Society 64, 477-480 (1942)

2. Lee HM., Wood DM., Hudson S., Archer JRH., Dargan PI. *Journal Medicine Toxicology* 10, 299-302. (2014)
3. Welter J., Meyer MR., Wolf EU., Weinmann W., Kavanagh P., Maurer HH. *Analytical Bioanalytical Chemistry* 405, 3125-3135 (2013)
4. WADA Technical Document – TD2017MRPL, Minimum required performance levels for detection and identification of non-threshold substances (2017)
5. WADA Technical Document TD2015IDCR. Minimum Criteria for Chromatographic-Mass Spectrometric Confirmation of the Identity of Analytes for Doping Control Purposes. (2014).

WP-107 / QUANTITATION OF FIVE BARBITURATES IN HUMAN URINE BY LIQUID CHROMATOGRAPHY-TANDEM MASS-SPECTROMETRY

Anna Chulovskaya - Dmitry Rakov - Elena Mochalova

National Anti-Doping Laboratory, MSU, Moscow

Keywords: barbiturates, validation, UPLC-MS/MS, forensic toxicology

Introduction

Barbiturates are a group of medicinal substances with a wide spectrum of central nervous system inhibitory effects, from mild sedation to total anesthesia [1]. They have a similar chemical structure to each other. Here, we present rapid, accurate and sensitive UPLC-MS/MS methods for the simultaneous quantitation of five barbiturates (butobarbital, butalbital, pentobarbital, amobarbital and secobarbital) in human urine.

Methods

Compounds were extracted from human urine samples using the LLE. The sample analysis was performed using Thermo Ultimate 3000 UPLC coupled with TSQ Vantage triple quadrupole mass spectrometer. Chromatographic separation was performed on an Acquity UPLC BEH C18 (2.1 mm × 100 mm, 1.7 μm) column coupled with Acquity UPLC BEH C18 (2.1 mm × 50 mm, 1.7 μm) guard-column. The 0.1 % formic acid in water and acetonitrile were used as A and B mobile phase, respectively.

Results

The main difficulty behind this study is to develop a method, allowing to distinguish structural isomers pentobarbital and amobarbital. The appropriate compounds separation was achieved by linear gradient from 27 % to 95 % for 8.5 min in 0.35 ml/min at 60°C. The method was fully validated for human urine in the range of 50 – 1000 ng/ml. Correlation coefficients (r^2) for five various barbiturates were > 0.99. Recovery, repeatability, and within-laboratory reproducibility lie in the range of 89-97%, 3-5%, and 7-8%, respectively.

Conclusions

A rapid and robust UPLC-MS/MS method has been developed and validated for the simultaneous determination of butobarbital, butalbital, pentobarbital, amobarbital and secobarbital in human urine. We proposed the LC conditions for the pentobarbital and amobarbital separation which allows to differ these compounds. The method are suitable for quantitation of barbiturates not only in toxicology, but in pharmacokinetic studies and clinical trails.

Novel Aspect

High chromatographic separation power for the specificity of five various barbiturates by UPLC-MS/MS in human urine was achieved in the field of forensic toxicology.

References

1. Zhang X., Lin Z., Li J., Huang Z., Rao Y., Liang H., Yand J., Zhenga F., *Drug Test Anal.* 9, 588-595 (2015)

WP-108 / A SYSTEMS BIOLOGY APPROACH FOR DISCOVERING BIOMARKERS INDICATIVE OF EXPOSURE TO CARFENTANIL

Elizabeth Dhummakupt⁽¹⁾ - *Bao Tran*⁽¹⁾ - *Daniel Carmany*⁽²⁾ - *Gabrielle Boyd*⁽²⁾ - *Phillip Mach*⁽²⁾ - *Michael Feasel*⁽¹⁾ - *Jennifer Sekowski*⁽¹⁾ - *Trevor Glaros*⁽¹⁾

Edgewood Chemical Biological Center, Research & Technology Directorate, Aberdeen Proving Grounds (1) - Excet, Inc., Research & Technology Directorate, Springfield (2)

Keywords: opioids, metabolomics, proteomics, biomarkers

Introduction

There are very few validated assays that are field-capable of rapidly detecting carfentanil, an ultra-potent, synthetic opioid related to fentanyl. Carfentanil has been responsible for numerous reports of mass localized overdoses in the United States. In order to improve the ability to detect this compound, a systems biology approach was utilized to discover host-based markers that are amplified upon carfentanil exposure.

Methods

New Zealand white rabbits were exposed to a 0.1 ECT50 concentration of carfentanil through nose-only inhalation. Blood was drawn immediately before and after exposure, then at time points out to thirteen days post-exposure. Blood plasma was then analyzed utilizing proteomics and untargeted metabolomics pipelines. For proteomics, a multidimensional fractionation method was developed for deeper coverage, and isobaric tandem mass tags (TMT) were utilized for quantitative assessments.

Results

For the proteomics, each animal served as its own control, due to the high inter-animal variability. By utilizing this method and stringent data processing filters, over 700 proteins were selected for statistical analysis. Of those, over 40 proteins were identified to have a two-fold ratio change in at least one time point, and 19 proteins resulted in significant p-values from ANOVA analysis. The metabolomics analysis identified over 100 features that had both a two-fold ratio change and a significant p-value from ANOVA. The PLS-DA showed significant separation when comparing pre- and post-exposure to the six-hour time point. Separation continued to become more pronounced as time progressed in both PLS-DA and PCA; however, after the first day the separation appeared to be maximized out to the latest time point. Identified proteins and metabolites suggest carfentanil exposure causes increased reactive oxygen species, possible vascular damage from coagulation and clotting, and possible myocardial infarction.

Conclusions

Proteomics and metabolomics pipelines were leveraged for a long time course study of low-dose exposure to carfentanil. A set of 19 proteins and over 100 metabolite features were identified as significant via fold-change and ANOVA. Both sets of data indicated significant changes could be seen six to 24 hours after exposure. Remarkably, a handful of metabolic perturbations were observed in samples taken immediately post exposure.

Novel Aspect

First utilization of proteomic and metabolomic pipelines for time course analysis following carfentanil exposure.

Funding for this project is provided by the Defense Threat Reduction Agency – Joint Science and Technology Office for Chemical and Biological Defense.

WP-109 / CO-EXPOSURE TO SILVER NANOPARTICLES AND CADMIUM IONS INDUCES METABOLIC ADAPTATION IN HEPG2 CELLS

Renata Miranda⁽¹⁾ - **Vladimir Gorshkov**⁽¹⁾ - **Barbara Korzeniowska**⁽¹⁾ - **Stefan Kempf**⁽¹⁾ - **Francisco Neto**⁽²⁾ - **Frank Kjeldsen**⁽¹⁾

University of Southern Denmark, Odense (1) - Federal University of Parana, Curitiba (2)

Keywords: Nanotoxicology, Proteomics, Silver nanoparticles, In vitro liver cell culture

Introduction

Although studies have reported the toxicological effects and underlying mechanisms of toxicity of silver nanoparticles (AgNP) in a variety of organisms, the interactions of AgNP with environmental contaminants are poorly understood. We used biochemical assays and mass spectrometry-based proteomics to assess the molecular effects induced by a co-exposure of HepG2 cells to AgNP and cadmium ions.

Methods

Human HepG2 hepatocellular carcinoma cells (Sigma-Aldrich) were cultured as a monolayer in high-glucose Dulbecco's modified Eagle's medium (DMEM) and subjected to spherical AgNP (10 nm), cadmium, or a combination of the two for 4-h and 24-h's for comparison. Cell viability, ATP/ADP ratio assays were performed in parallel to quantitative proteomics (Q Exactive HF) using isotopic labeling.

Results

Cell viability and energy homeostasis were slightly affected after 4-h exposure to AgNP, Cd²⁺, or a combination of the two; these endpoints were substantially altered after a 24-h co-exposure to AgNP and Cd²⁺, while exposure to one of the two contaminants led to minor changes. Proteomics analysis followed the same trend: while a 4-h exposure induced minor protein deregulation, a 24-h exposure to a combination of AgNP and Cd²⁺ deregulated 43% of the proteome. The toxicity induced by a combined exposure to AgNP and Cd²⁺ involved (1) inactivation of Nrf2, resulting in downregulation of antioxidant defense and proteasome-related proteins, (2) metabolic adaptation and ADP/ATP imbalance, and (3) increased protein synthesis intended to reestablish homeostasis. The adaptation strategy was not sufficient to restore ADP/ATP homeostasis and to avoid significant cell death.

Conclusions

Co-exposure for 24 h to AgNP and Cd²⁺ induced more toxic responses than individual exposures. The toxicity induced by AgNP+Cd²⁺ was the upregulation of oxidative phosphorylation and lipid metabolism, suggesting a metabolic adaptation, as HepG2 cells usually produce ATP mainly by glycolysis, and alternative nutrients, such as lipids, may have been recruited to reestablish an energy homeostasis.

Novel Aspect

Our results provide information on the toxicological interaction of AgNP and Cd²⁺ ions in HepG2 cells, and insights into possible toxicity mechanisms.

WP-110 / APPLICATION OF DISPERSIVE LIQUID-LIQUID MICROEXTRACTION (DLLME) FOR IDENTIFICATION OF SYNTHETIC CANNABINOIDS IN PLASMA BY LC-MS/MS

Lino Thais⁽¹⁾ - **Ingrid Lopes Barbosa**⁽²⁾ - **Karina Diniz de Oliveira**⁽³⁾ - **Marcos N. Eberlin**⁽²⁾ - **José Luiz Costa**⁽⁴⁾

School of Medical Sciences, University of Campinas, Rua Tessalia Vieira de Camargo, 126, Campinas, São Paulo, SP, Brazil (1) - Thomson Mass Spectrometry Laboratory, Chemistry Institute, University of Campinas, R. Sérgio Buarque de Holanda, 126 - Cidade Universitária, Campinas, São Paulo, SP, Brazil (2) - Division of

Psychiatry, School of Medical Sciences, University of Campinas, Rua Tessalia Vieira de Camargo, 126, Campinas, São Paulo, SP, Brazil (3) - Faculty of Pharmaceutical Sciences, University of Campinas, Campinas (4)

Keywords: Synthetic cannabinoids, DLLME, LC-MS/MS, New Psychoactive Substances

Introduction

The detection of synthetic cannabinoids (SCs) in biological fluids is challenging and require sensitive analytical methods. The dispersive liquid-liquid microextraction (DLLME) is a sample preparation technique based on the partition of analytes using small volumes of solvent mixture (dispersive and extractor). The aim of this work is optimize a DLLME extraction to analyze SCs in plasma by LC-MS/MS.

Methods

Sample (400 µL) was fortified with 40 µL of SCs mixture, vortexed and 1200 µL of dispersive solvent was added and centrifuged at 5500 rpm/5 min. Water (1 mL) and sodium tetraborate saturated solution (100 µL) was added to 500 µL of supernatant and 100 µL of the extraction solvent was rapidly added and the mixture was centrifuged at 5500 rpm/5 min. An aliquot (180 µL) was transferred to a vial, evaporated, reconstituted in 100 µL of methanol and injected (2 µL) into LC-MS/MS.

Results

The dispersive and extraction solvents were rationally chosen by design of experiment (D-optimal), based on their solubility and extraction capacity, capacity to the formation of the cloudy (dispersive effect). The dispersive solvent has to be miscible with both aqueous and organic phases, in order to ensure the formation of the cloudy state that enhances the contact between two phases, thus facilitating extraction. The design of experiment pointed out that the best solvents for extraction was acetonitrile as dispersive solvent, and chloroform as extractor solvent. This combination of solvents allowed the efficient extraction of SCs as low as 0.1 ng/mL, using small plasma volume (0.4 mL). DLLME still needs to be improved and further experiments will be carried out, to evaluate the optimum condition of sample volume, pH, salting-out process and centrifugation conditions.

Conclusions

The extraction method initially optimized and applied for the extraction of SCs was adequate to recover the analytes in plasma samples at low concentrations (0.1 ng/mL), concentration expected to be found in real cases. We're improving our method, performing other studies to improve even more the extraction steps, and thus the method detectability, to finally validate and apply in real samples.

Novel Aspect

DLLME is a microextraction technique, which requires low volume of solvents, and was successfully applied to SCs analysis, extremely toxic emergent drugs.

WP-111 / DEVELOPMENT AND VALIDATION OF AN LC-MS/MS METHOD FOR THE DETERMINATION OF MATERNAL AND FETAL TISSUE DISTRIBUTION OF PERMETHRIN AND ITS METABOLITES IN PREGNANT RATS

Stephane Personne⁽¹⁾ – *Aurélie Corona*⁽²⁾ - *Celine Brochot*⁽³⁾ - *Karen Chardon*⁽¹⁾ - *Sylvie Baltora*⁽⁴⁾ - *Serge Pilard*⁽⁵⁾ - *Veronique Bach*⁽¹⁾ - *Florence Zeman*⁽³⁾ - *Paulo Marcelo*⁽⁶⁾

Universite de Picardie Jules Verne, UMRI 01 INERIS PERITOX, Amiens (1) - INERIS, Unite TOXI, Verneuil en Halatte (2) - INERIS, Unite METO, Verneuil en Halatte (3) - Universite de Picardie Jules Verne, UMR7058 CNRS, Amiens (4) - Universite de Picardie Jules Verne, Plateforme PFA, Amiens (5) - Universite de Picardie Jules Verne, Plateforme ICAP, Amiens (6)

Keywords: UHPLC-MS/MS, Permethrin, Toxicokinetic, Gestation, Rat

Introduction

Animal studies suggest that gestational exposure to permethrin can cause neurodevelopmental toxicity [1]. Permethrin can cross the placental barrier [2]. However, toxicokinetic data are not available to quantify this transfer.

The aim of our study was to develop an UHPLC-MS/MS method to quantify permethrin isomers and their metabolites in blood, urines, feces and several tissues of interest: brain, liver, kidney, mammary gland, fat and placenta in pregnant rats.

Methods

The extraction was performed by liquid/liquid extraction. Chromatographic conditions for cis and trans-permethrin and their metabolites (cis and trans-DCCA and 3-PBA) were optimized using different analytical columns, mobile phases and way of elution.

The final method was validated for recovery, specificity, linearity, precision and accuracy before its application to quantify samples of a toxicokinetic study in pregnant rats exposed daily to permethrin.

Results

Chromatographic separation was accomplished on a HSS T3 column with a gradient elution system. Two different mobile phases were used. For permethrin, the mobile phase consisted of ammonium acetate buffer and acetonitrile. For metabolites, acetonitrile and water containing 0.1% of formic acid were selected.

Relative recoveries were found to be in the range of 66 to 117% for permethrin and 43% to 140% for metabolites. No interference from other components in the matrix were observed at the retention times of the compounds.

Matrix-matched calibration curves had satisfactory linearity up to 2000 ng/mL. The LOQ ranged from 10 to 200 ng/mL for the compounds depending on the matrix. The intra- and inter-batch precision and accuracy were better than 15%.

The validated method was applied to analyze samples collected in the toxicokinetic study. Cis-permethrin, trans-DCCA and 3-PBA were detected in all matrices, including fetal matrices. Trans-permethrin and cis-DCCA were detected in all matrices except placenta and liver respectively.

Conclusions

We developed a rapid analytical method with limited sample pre-treatment sufficient to allow the simultaneous quantification of permethrin and their metabolites in biological matrices in pregnant rats and fetuses.

The method showed good performances in terms of accuracy and precision. The LOQ were better compared to a previously published method targeting the same compounds in a large panel of biological matrices in rats [3].

Novel Aspect

This method is the first method used to successfully quantify permethrin isomers and their metabolites in a various panel of biological matrices in pregnant rats and fetuses.

References

1. Shafer TJ., Meyer DA., Crofton KM. Critical Review and Future Environ Health Perspect, 113, 123-136 (2004).
2. Pérez JJ., Williams MK., Weerasekera G. J Chromatogr B, 878, 2554-2562 (2010).
3. Lestremau F., Willemin M-E., Chatellier C., Desmots S., Brochot C. Anal. Bioanal. Chem. 406, 3477-3487 (2014).

WP-112 / PROTON TRANSFER REACTION – MASS SPECTROMETRY AND ITS APPLICATIONS TO HOMELAND SECURITY: DETECTION OF ILLICIT DRUGS

David Olivenza-Leon⁽¹⁾ - *Peter Watts*⁽¹⁾ - *Fraser Reich*⁽²⁾ - *Chris Mayhew*⁽³⁾

University of Birmingham, School of Physics and Astronomy, Birmingham (1) - Kore Technology Ltd, Kore Technology Ltd, Ely (2) - Universitat Innsbruck, Institut Fur Atemgasanalytik, Innsbruck (3)

Keywords: Proton transfer reaction mass spectrometry, illicit drugs, thermal desorption, homeland security

Introduction

Proton Transfer Reaction Mass Spectrometry (PTR-MS) is an analytical technique for the detection of trace chemical compounds in complex environments through soft chemical ionization useful in many fields including health sciences, environmental sciences and homeland security. PTR-MS is based on the volatility of the targeted chemicals and its ability to work without initial separation or preconcentration makes it suitable for on-line, real time measurements.

Methods

PTR-MS is a suitable analytical technique to monitor volatile organic compounds (VOCs) in the air. However, compounds of interest in the field of homeland security can have low pressure, which makes their detection challenging. To amend this, a Thermal Desorption Unit (TDU) in conjunction with a PTR-MS instrument can be used. It has recently been applied to the detection of explosives [1] and here we apply it to the detection of illicit drugs.

Results

In this study we demonstrate the application of a TDU for the on-line detection of trace quantities of the most widely used societal drugs, the reactions of which with H_3O^+ at just one reduced electric field, 120 Td, have already been reported [2].

The results shown were achieved by measuring some hundreds of nanograms deposited in a teflon swab which was subsequently inserted into the TDU for the desorption of the sample into the PTR-MS instrument. The two modes of the instrument, namely DC and RF modes, were used. In DC mode, a uniform electric field draws the ions downstream in the reactor, while in RF mode a radio frequency field is superimposed to the uniform electric field to deliver extra collisional energy to the molecules.

Conclusions

This study shows that the TDU is a suitable tool to use in conjunction with PTR-MS for the detection of illicit drugs. The relative distribution of product ions in the DC mode agrees with the literature [2], but the results from the RF mode is controversial. For heroin, morphine and codeine both modes present similarities, showing more fragmentation of the parent ion in the RF mode, whereas for cocaine and MDMA they totally differ.

Novel Aspect

The novelty of these results lies in the use of the TDU and the RF mode in a PTR-MS instrument for the detection of illicit drugs.

References

1. González-Méndez R., Reich D., Mullock S., Corlett C., Mayhew C., International Journal of Mass Spectrometry 385, 13-18 (2015)
2. Agarwal B., Petersson F., Jürschik S., Sulzer P., Jordan A., Märk T., Watts P., Mayhew C., Analytical and Bioanalytical Chemistry 400, 2631-2639 (2011)

Aleksandra Pawlaczyk - Martyna Balcerek - Małgorzata Iwona Szyrkowska

Lodz University of Technology, Faculty of Chemistry, Institute of General and Ecological Chemistry, Lodz

Keywords: trace elemental analysis, eye shadows, cosmetics, ICP-MS, CV-AAS

Introduction

Among color cosmetics eyeshadows deserve special attention since they are applied in the area where facial skin is thin. Many studies suggest that cosmetics belonging to the eye makeup segment contain high amount of toxic metals. The size of particles used in their production also plays a key factor in safety issues. All may contribute to the increasement of allergic reactions risk or development of other sides effects appearing with their prolonged application [1-2].

Methods

ICP Mass Spectrometry with Time-of-Flight Analyzer (Ag, Ba, Bi, Cd, Pb, Sr and Tl) and Cold Vapor Atomic Absorption Spectrometry (Hg) techniques were applied for inorganic chemical elements determination in eyeshadows products. In total 94 samples (in 12 colors) obtained from the local market and originating from 8 most commonly used brands were purchased for the analysis. The samples were coded depending on their price, color, type, origin and producer.

Results

They are no current data concerning the potential assessment of the risk connected with application of eyeshadows. The safety of their use is dependent upon *e.g.* the type of products (source and form of the components, size of the particles), time and frequency of the application, etc. In this study eyeshadows samples were grouped by similarity to the concentrations of the chosen chemical elements obtained by ICP-MS and CVAAS techniques. The price of products, their color, brand, origin, destiny (for kids and adults), type (matt and pearl) were used as criteria of their classification. Multivariate methods as PCA and CA were applied to observe the relationships and interactions of several analyzed elements and assess influence of chosen factors. It was shown that the average level of Cd in one sample was few times above the permissible limit (0.5 mg/kg), while the concentration of Pb and Hg in two samples should not be ignored since they slightly exceeded the polish regulations regarding cosmetics products (10 mg/kg and 1 mg/kg, respectively).

Conclusions

Summarizing the influence of chosen factors it can be concluded that factors having the least impact on the concentration of studied elements are color and type of eyeshadows. Factors affecting the majority of concentrations of the determined elements are the brand of the eyeshadow (Sr, Ag, Cd, Pb, Bi, Hg), country of origin (Sr, Ba, Pb, Bi, Hg) and shadow price (Sr, Ba, Pb, Bi, Hg). For the following elements: Sr, Pb, Bi, Hg the influence of the same factors was observed.

Novel Aspect

The examination of the content of chosen trace metals in the locally sourced eyeshadows is quite relevant since there are no updated reports regarding homegrown products safety.

References

1. Santos B. V., Martins de Oliveira Junior J., Bonventi Junior W., Miwa Hanai-Yoshida V., X-Ray Spectrometry 47, 242 (2018).
2. Contado C., Pagnoni A., Science of the Total Environment 432, 173 (2012).

WP-114 / GENERIC EXTRACTION METHOD FOR SCREENING OF 8 NARCOTIC DRUGS IN SALIVA WITH A QUICK TURNAROUND ANALYSIS BY LDTD-MS/MS

Pierre Picard - Jean Lacoursiere - Serge Auger

Phytronix Technologies Inc, 4535 boul. Wilfrid-Hamel, suite 120, Quebec

Keywords: Narcotics, Saliva, LDTD-MS/MS, screening

Introduction

A generic extraction method combined with LDTD-MS/MS analysis was developed for fast turnaround screening of drugs in saliva. This new method could give police officers rapid and accurate answer in less than 10 minutes allowing on-site screening on a police roadblock.

Methods

Drugs are spiked at 50%, 100% and 150% of the decision point in saliva. 50 µL is mixed with a solution composed of 100 µL of internal standard and 75 µL of buffer. The mixture is vortexed 20 seconds and centrifuged 30 seconds. 8 µL of the upper layer are deposited in LazWell-HDE plates. The plates are dried for 3 minutes. Seven drugs are analyzed in positive APCI mode and one in negative APCI mode in a single run. LDTD's laser power was ramped from 0 to 55% of maximum power in 6 seconds.

Results

The decision points in saliva are 15 ng/mL for THC, 50 ng/mL for Amphetamine/Methamphetamine/MDMA and 10 ng/mL for Cocaine/Benzoylcegonine/Morphine/6-AM. Using the SWGTOX validation guidelines, a blank preparation of saliva is spiked with 8 previously enumerated drugs at concentration representing 50%, 100% and 150% of the cut-off. Validation of this new screening approach, the CV at each concentration must not exceed 20% and the average concentration ± 2 times the standard deviation must not overlap with other concentration at the decision point. Triplicate extractions are performed for each concentration over the course of 5 days, for validation. Additional blank of saliva, one with Internal standards (IS) and the others without are extracted and analyzed. The THC is analyzed in negative mode and the others in positive. No overlapping at the decision points is observed and a %CV lower than 15% is obtained. Series of tests were timed from the saliva collection to the output result giving an average of 9.3 minutes for the complete workflow.

Conclusions

Simple extraction method combined with LDTD-MS/MS analysis allow specific measurement of narcotic in saliva. The generic method enables accurate on-site analysis of saliva collected from an individual potentially driving under the influence of drugs with expected feedback results obtained under 10 minutes by the police officer.

Novel Aspect

In situ High Throughput targeted screening of narcotic in saliva using LDTD-MS/MS

WP-115 / IMPROVEMENT OF THE ATHLETE BIOLOGICAL PASSPORT BY BLOOD STEROID PROFILING: CHALLENGES OF A STEROIDOMIC APPROACH

Federico Ponzetto⁽¹⁾ - *Julien Boccard*⁽²⁾ - *Raul Nicoli*⁽¹⁾ - *Tiia Kuuranne*⁽¹⁾ - *Serge Rudaz*⁽²⁾ - *Martial Saugy*⁽³⁾

Swiss Laboratory for Doping Analysis, CHUV - University of Lausanne, Lausanne (1) - School of Pharmaceutical Sciences, University of Geneva, Geneva (2) - Center for Research and Expertise in anti-Doping sciences, University of Lausanne, Lausanne (3)

Keywords: steroidomics, doping, HRMS, serum

Introduction

The steroidal Athlete Biological Passport aims at revealing use of testosterone. It is implemented by GC-MS quantification of selected endogenous androgens in urine. Sensitivity of this approach suffers from confounding factors *e.g.* enzyme polymorphism, bacterial contamination and ethanol consumption [1]. Among the improvement strategies, quantification of steroid hormones in serum has shown promising results and a complement to the urinary module [2, 3].

Methods

A steroidomic approach was used to highlight new serum biomarkers of testosterone doping. The developed analytical workflow consisted of a SPE clean-up of steroids and their phase II metabolites, followed by UHPLC-HRMS analysis. Full Scan MS and data dependent MS/MS spectra were acquired in both positive and negative ionization modes. Samples from testosterone administration study (n=380) were analyzed in 11 analytical batches.

Results

Positive and negative ion mode data from UHPLC-HRMS analyses were treated separately: firstly, data processing including noise filtering, chromatographic alignment and peak picking was performed by means of dedicated software and resulted in datasets with approximately 40000 features for each ionization mode. Then, in order to reduce data dimensionality, ion annotation was carried out assigning confidence level 1 (confirmed with authentic standard) and 2 (consistent with predicted properties) employing a newly created in-house steroid database (350 steroids) together with retention time prediction model [4]. Intensities of annotated features (64 and 174 in positive and negative mode, respectively) were corrected with LOESS regression in order to minimize the batch effects and finally Analysis of variance Multiblock Orthogonal Partial Least Squares (AMOPLS) was performed to highlight new promising blood biomarkers of testosterone intake, similarly to what was previously done on urine matrix [5].

Conclusions

The study revealed a set of biomarkers for both oral and transdermal testosterone abuse, including hydroxylated metabolites of androgens, as well as their glucuronidated phase II metabolites. In comparison with urinary steroid profiling, these targets could increase sensitivity and extend the current detection window. In the near future, targeted quantitative UHPLC-MS/MS method will be developed for serum to validate the results of this study.

Novel Aspect

This work represents a first application of steroidomics on serum samples for anti-doping purposes, demonstrating the usefulness of holistic approaches in such a context.

References

1. Kuuranne T., Saugy M., Baume N., Br. J. Sports. Med., 48, 848-855 (2014).
2. Thevis M., Kuuranne T., Walpurgis K., Geyer H., Schänzer W., Drug Test. Anal., 8, 7-29 (2016).
3. Ponzetto F., Mehl F., Boccard J., Baume N., Rudaz S., Saugy M., Nicoli R., Anal. Bioanal. Chem., 408, 705-719 (2016).
4. Randazzo G. M., Tonoli D., Hambye S., Guillarme D., Jeanneret F., Nurisso A., Goracci L., Boccard J., Rudaz S., Anal. Chim. Acta, 916, 8-16 (2016).
5. Boccard J., Badoud F., Jan N., Nicoli R., Schweizer C., Pralong F., Veuthey J. L., Baume N., Rudaz S., Saugy M., Bioanalysis, 6, 2523-2536 (2014).

WP-116 / ANALYSIS OF REMOVAL OF PESTICIDES AND DRUGS IN WASTE WATER BY PLANTS USING FLOW INJECTION ANALYSIS (FIA) MAGNETIC RESONANCE MASS SPECTROMETRY (MRMS)

***Matthias Witt*⁽¹⁾ - *Aiko Barsch*⁽²⁾ - *Claire Villette*⁽³⁾ - *Dimitri Heintz*⁽³⁾**

Bruker Daltonik GmbH, MRMS Solutions, Bremen (1) - Bruker Daltonik GmbH, Metabolomics, Bremen (2) - IBMP-CNRS, University of Strasbourg, Strasbourg (3)

Keywords: Pesticides, Flow Injection Analysis, Drug detection, Screening

Introduction

Treatment of wastewater is important to remove pollutants, drugs and other potentially environmentally toxic molecules before release into natural compartments. Constructed wetlands in France are used to treat water from small villages. Here we present the profiling of plants to understand if pesticides and drugs can be removed from waste water through an accumulation in plants or degradation by environmental factors.

Methods

Two poplars (*Populus nigra*) were planted - one close to the border of a pond and the other away from the pond (control). Mature leaves of these two plants were collected and extracted. Each sample was analyzed in 8 technical and 3 measurement replicates in FIA-ESI using positive and negative ion mode. A solarix 2xR (Bruker Daltonik GmbH, Bremen) MRMS system was used for these studies. Analysis of data was performed with MetaboScape 3.0 (Bruker Daltonik GmbH, Bremen).

Results

The data of the ESI(+) and ESI(-) were combined for feature analysis. More than 3400 features have been found for the plant extract samples. Roughly 90% of the detected features could be assigned with a molecular formula using a mass tolerance of only 0.5 ppm. More than 360 compounds have been annotated in the samples using a food data base with nearly 16000 entries using same mass tolerance. Features responsible to differentiate treated and non-treated plant sample were detected by T-test and Principle Component Analysis (PCA). These compounds were analyzed in detail as possible drugs and pesticides. Based on these results 6 pesticide and 4 drug candidates were found by screening versus pesticide and drug data bases.

Conclusions

Pesticides and their metabolites can be detected by FIA-MRMS in plant samples. Statistical analysis such as T-test and PCA can be used to find compounds responsible to distinguish treated and non-treated plants based on FIA-MRMS data. This method can be used to understand if pesticides and drugs can be removed from waste water.

Novel Aspect

FIA-MRMS can be used to study the effect of removal of drugs and pesticides in waste water.

WP-117 / IMPROVEMENT FOR HIGH SENSITIVITY OF THE DRUG SCREENING BY THERMAL DESORPTION AND PYROLYSIS COMBINED WITH DART-MS (TDP/DART-MS)

Hiroko Abe⁽¹⁾ - *Chikako Takei*⁽²⁾ - *Kenichi Yoshizawa*⁽²⁾ - *Motoshi Sakakura*⁽³⁾ - *Teruhisa Shiota*⁽³⁾ - *Daisuke Yajima*⁽¹⁾ - *Hiroto Iwase*⁽⁴⁾

University of Chiba, Department of Legal Medicine, Chiba (1) - BioChromato, Inc., Application support, Fujisawa, Japan (2) - AMR Inc., Application, Meguro-ku (3) - University of Tokyo, Department of Forensic Medicine, Bunkyo-ku (4)

Keywords: drug screening, urine, TDP/DART-MS

Introduction

In recent years, new drugs appear one after another, including new psychoactive substance that threaten

society. Comprehensive analytical method for drug screening is greatly desired, since it is difficult to examine individually the analytical conditions that are appropriate for each one. Thus, we are developing a comprehensive drug screening method for urine or blood sample without pretreatment by thermal desorption and pyrolysis-DART-MS (TDP/DART-MS).

Methods

Standard drug mixture solution and urine added drugs (*i.e.*, blank urine with several kinds of drug mixture added) were used as samples. TDP/DART-MS were composed by the combination TripleTOF 5600+ (Sciex) Q-TOF mass spectrometer equipped DART ion source (IonSense) and ionRocket TDP device (BioChromato). Mass spectra were measured as samples were heated from room temperature to 300°C at a rate of 60°C/min, the measurement time was 5 min.

Results

At first, analysis of standard drug mixture by TDP/DART-MS, each drug was separated and detected through adding gradient heating, and thermal desorption profiles were highly reproducible for individual drug. Additionally, the detected ions were correctly identified according to their measured accurate mass and MS/MS spectra. So, it was confirmed that TDP/DART-MS can be a suitable analytical method for our investigation. Next, we examined the minimum pretreatment such as deproteinization by ACN, since for analysis of urine added drugs without any pretreatment, it was indicated that we had to improve the detection sensitivity. Comparing non-pretreated and deproteinization pretreated, the detection intensity of each drug in urine sample pretreated deproteinization were increased considerably. Moreover, the calibration curves using urine were linear in that ranging 0.01–1 µg/ml.

Conclusions

In summary, TDP/DART-MS can identify drugs in urine, directly and rapidly. In addition, since the calibration curves were linear in that ranging 0.01–1 µg/ml. Although our issue is to improve the detection intensity of drugs, TDP/DART-MS will be useful for rapidly quantitative analysis of drugs in urine in the future, and this analysis method will be a useful method in the field of required rapidity such as hospital and police.

Novel Aspect

TDP/DART-MS is useful analysis method for rapid identification and screening of drugs in urine samples.

WP-118 / DETECTION OF IGF-1 IN SERUM AS A BIOMARKER OF GROWTH HORMONE ABUSE BY ATHLETES

Maxim Yunin - Irina Zvereva - Dmitry Rakov - Elena Mochalova

National Anti-Doping Laboratory, MSU, Moscow

Keywords: IGF-1, trypsin digestion, LC-MS/MS, doping-control

Introduction

The measurement of concentrations of insulin-like growth factor-1 (IGF-1) and *N*-terminal pro-peptide of type III collagen in serum was suggested for detection of growth hormone abuse by athletes [1]. The method for IGF-1 concentration determination by the quantification of two “signature peptides” (T1, T2) by UHPLC-MS/MS after trypsinolysis was proposed [2]. Here, we describe our results obtained after optimization and validation of the above method.

Methods

T1 and T2 were obtained after IGF-1 binding protein complex dissociation followed by protein precipitation by acetonitrile and trypsinolysis. Calibrations curves were made using rat serum fortified with human IGF-1. Also we used calibrators that were stored at -80°C. T1 and T2 peptides were analyzed by Thermo Ultimate

3000 UHPLC with TSQ Vantage triple quadrupole mass spectrometer. The structures of peptides were confirmed by HRMS on Q Exactive mass spectrometer.

Results

The range of bottlenecks was revealed during method implementation. Calibration solutions that allow the estimation of the inter-day precision show remarkable changes in the standard concentration over time. The critical step in sample preparation was DTT reduction. Incubation for more than 90 min at 60°C led to T1 signal disappearance. Thus, we suppose that milder DTT reduction conditions are apparently more preferable for working with long peptides. For UHPLC-MS/MS method, the use of the gentle gradient (5-45%B, 20 min, 0.22 mL/min) at 40°C and high dwell time (200 ms) were necessary for achieving a good response at low levels. The method was fully validated in the range of 50 – 600 ng/mL. The correlation coefficients (r^2) for T1 and T2 were $0.996\pm 0.004/0.996\pm 0.003$, respectively. Interestingly, the correlation coefficients (r^2) for T1 and T2 in the range of 50-1000 ng/mL were $0.984\pm 0.003/0.984\pm 0.003$, respectively. Therefore, the calibration point at 1000 ng/mL of IGF-1 was recognized as outlier and removed.

Conclusions

The method was fully validated and implemented for anti-doping purposes. Determination of T1 peptide was the critical point of the method, probably, due to its longer amino acid sequence as compared to T2. It was determined that high temperatures should be avoided during the reduction step. For LC conditions, the gentle gradient, low flow rate and column temperature should be used to facilitate separation T1. High dwell time is crucial for method sensitivity.

Novel Aspect

The incubation conditions in the presence of DTT, as well as UHPLC-MS/MS parameters have a significant effect on the T1 peptide sensitivity.

References

1. World Anti-Doping Program: Human growth hormone (hGH) biomarkers test for doping control analysis (2016), April 2016 (version 2).
2. Cox H. D., Hughes C. M., Eichner D., *Bioanalysis*, 6(19), 2651-2662 (2014).

WP-119 / DEVELOPMENT OF A LC-HRMS BASED STRATEGY FOR THE SCREENING OF CONJUGATED STEROIDS IN EQUINE URINE

Xuan Zuo⁽¹⁾ - **Zied Kaabia**⁽¹⁾ - **Florian Bagilet**⁽¹⁾ - **Béatrice Georges**⁽¹⁾ - **Marie Collard**⁽¹⁾ - **Yves Moulard**⁽²⁾ - **Patrice Garcia**⁽³⁾ - **Marie-Agnès Popot**⁽³⁾ - **Ludovic Bailly-Chouriberry**⁽²⁾

Laboratoire des Courses Hippiques (LCH), Anabolic agents anti-doping department, Paris, France (1) - Laboratoire des Courses Hippiques (LCH), Medication control department, Paris, France (2) - Laboratoire des Courses Hippiques (LCH), Research and Development department, Paris, France (3)

Keywords: LC-HRMS, Androgenic anabolic steroids, Horse racing, equine urine, Q-Exactive HF-X

Introduction

The gold standard method for analyzing androgenic anabolic steroids (AAS) in equine urine, in the context of anti-doping control, is GC-MS/MS after deconjugation and derivatization steps. LC-HRMS can offer an alternative option for the direct analysis of AAS under their intact conjugated forms in horse urine. The aim of this study is to develop an analytical strategy allowing the extraction of conjugated AAS and their analysis by LC-HRMS using the Q-Exactive HF-X.

Methods

AAS are mainly under their sulfo or glucuruno-conjugated forms in equine urine. Therefore an adapted

analytical strategy was settled for their detection. The sample preparation was optimized assessing different SPE cartridges and conditions, which have been compared in terms of extraction recovery. The chromatographic separation and the mass spectrometry parameters were also optimized for conjugated steroids detection with higher sensitivity and selectivity.

Results

SPE cartridges with C18 or strong anion exchange (SAX) stationary solid phases were evaluated during the study. Extraction of conjugated AAS from horse urine were performed on C18 cartridges using percentages of methanol in wash steps ranging between 5 and 40%. Higher percentages of methanol resulted in cleaner extracts and lower matrix noise impact. However, increasing this percentage induced a loss of some glucuronide and di-glucuronide steroids. Unlike the C18 cartridges, the anion exchanger cartridges resulted in a higher extraction recovery of glucuronide steroids and lower extraction yields for sulfo-conjugated forms. The selected chromatographic condition using C18 column showed well-resolved chromatographic peaks and a good separation profile. Finally, the MS parameters were optimized using the ESI negative mode to reduce background noise and enhance conjugated steroids ionization. A scan mode combining Full MS and PRM was settled and AGC parameters optimized to enhance sensitivity and selectivity.

Conclusions

The developed analytical strategy was based, first, on a steroid extraction performed using a C18 SPE cartridge characterized by an efficient extraction of conjugated steroids from horse urine together with an optimized washing step reducing matrix effect during the LC-MS analysis. Then, the optimization of chromatographic and mass spectrometry parameters permitted to obtain satisfying separation and detection of sulfo and glucuruno-conjugated forms.

Novel Aspect

Unlike gold standard GC-MS/MS strategy, LC-HRMS allows conjugated AAS analysis in equine urine under their intact form insuring a gain in time and avoiding potential artifacts.

References

1. International agreement on breeding, racing and wagering. International Federation of Horseracing Authorities (IFHA), Paris, May 2018.

WP-120 / FAST GC-MS/MS FOR HIGH THROUGHPUT ANDROGENIC ANABOLIC STEROIDS ANALYSIS IN EQUINE URINE

Karen Histel⁽¹⁾ - *Marie Collard*⁽¹⁾ - *Béatrice Georges*⁽¹⁾ - *Florian Bagilet*⁽¹⁾ - *Marie-Agnès Popot*⁽²⁾ - *Zied Kaabia*⁽¹⁾ - *Ludovic Bailly-Chouriberry*⁽³⁾

Laboratoire des Courses Hippiques (LCH), Anabolic agents anti-doping department, Paris, France (1) -
Laboratoire des Courses Hippiques (LCH), Research and Development department, Paris, France (2) -
Laboratoire des Courses Hippiques (LCH), Paris, France (3)

Keywords: Fast GC, Androgenic anabolic steroids, Horse racing, equine urine, mass spectrometry

Introduction

The use of androgenic anabolic steroids (AAS) as growth promoters for equine performance enhancement is forbidden according to article 6 of the International Federation of Horseracing Authorities (IFHA) [1]. The latter bans the use of substances capable of giving horses an advantage, contrary to the horse's inherent merit. Moreover, the use of these substances is supposed to result in biological and behavioral modifications, considered dangerous for the animal.

Methods

The gold standard method for analyzing AAS in equine urine is GC-MS/MS following hydrolysis and extraction steps. The duration of the analytical run is circa 30 minutes per sample which consists an issue for high throughput routine analysis. In order to reduce the analysis runtime and enhance performances, a Fast GC-MS/MS method has been developed using a shorter column of 20 m instead of 30 m with tighter i.d and thinner film in comparison to GC columns.

Results

Two Fast GC columns, Drug-1 (10 m x 0.18 mm x 0.18 μm) and ZB-5MS-PLUS (20 m x 0.18 mm x 0.18 μm), were evaluated and compared to the currently used GC column DB-5MS-UI (30 m x 0.25 mm x 0.25 μm). The chromatographic profile of the 47 steroids of the study, when performing the separation on the 10 m Drug-1 column, exhibited phenomena of peak tailing and carry-over. These phenomena were caused by the matricial charge density versus the low quantity of stationary phase. The use of a Fast GC column ZB-5MS-PLUS (20 m x 0.18 mm x 0.18 μm) with an optimized gradient permitted to separate the different steroids of the study and resulted in a well-resolved chromatographic profile. The Fast GC-MS/MS strategy using the 20 m Fast-GC column was performant in terms of repeatability, sensitivity and selectivity following the analysis of blank and spiked urine samples and reference standards of a 47 steroids pool mixture. The developed strategy led to a considerable gain in analysis time and hence, presents a significant interest for high throughput routine analysis.

Conclusions

The developed Fast GC-MS/MS method using the 20 m ZB-5MS-PLUS column allowed to obtain well-resolved chromatographic peaks for the different steroids of the study. The performance of the Fast-GC-MS/MS strategy were satisfactory namely in terms of sensitivity and robustness. In comparison to the currently used GC-MS/MS method, a gain in analysis runtime of 8 minutes per sample and 5 hours per instrument per day was obtained.

Novel Aspect

Fast-GC-MS/MS strategy results in a significant gain in analysis time and insures good performances which is crucial for high throughput routine analysis of AAS in equine urine.

References

1. International agreement on breeding, racing and wagering. International Federation of Horse racing Authorities (IFHA), Paris, May 2018.

WP-121 / IN-VITRO/VIVO STUDY OF ANDARINE METABOLISM, A SELECTIVE ANDROGEN RECEPTOR MODULATOR (SARM), IN THE CONTEXT OF HORSE DOPING CONTROL

Stéphane Trevisiol - Yves Moulard - Sophie Boyer - Murielle Jaubert - Benoit Loup - Patrice Garcia - Marie-Agnès Popot - Ludovic Bailly-Chouriberry

Laboratoire des Courses Hippiques, LCH, Verrière-le-buisson, France

Keywords: Doping, horse, metabolism, SARM, in-vivo/vitro

Introduction

Andarine is a selective androgen receptor modulator classified as non-steroid anabolic agent. Its use is forbidden by the International Federation of Horseracing Authorities at any time in horse career. This molecule have the ability to stimulate androgen receptors in muscles and bones and thus to enhance athletic performances. To control the fraudulent administration of andarine, a horse metabolism study has

been performed to identify the best target metabolites.

Methods

This metabolism study of andarine was carried out first by in-vitro approach using microsomal fractions of equine hepatocytes (S9) and then in-vivo on one horse via a subcutaneous administration route. Both in-vitro reaction products and blood/urine samples, collected up to 4 days, were analyzed by LC-HRMS2 and LC-SRM to characterize phase I and phase II andarine metabolites.

Results

The in-vitro results obtained using S9 equine hepatocyte fractions showed an important phase I metabolism of andarine. Indeed, six phase I andarine products were identified *i.e.* three metabolites with a functional loss (methyl, acetyl or phenyl), one with the nitro group reduction and one which was identified as cleavage product of andarine (4-nitro-3-trifluoromethyl-phenylamine) with its hydroxylated form. The post administration study has shown the presence of the native molecule in both urine and plasma. In urine, the results were consistent with these obtained in-vitro: the metabolite of andarine 4-nitro-3-trifluoromethyl-phenylamine and its hydroxylated form were also identified. These two metabolites were much longer detected in urine than the native molecule and presented high degree of phase II conjugation with both sulfate and glucuronic acid. Limit of detections (LOD) of andarine were 1 pg/mL in plasma and 5 pg/mL in urine. Concerning the two urine metabolites their LOD was around 50 pg/mL in andarine equivalent.

Conclusions

Andarine is prohibited at any time. Therefore, it was important to study its metabolism in horse to identify the best target metabolites to detect its misused. The in-vitro results demonstrated that andarine was strongly metabolized. These results have been confirmed by the in-vivo experiments. This study demonstrated that the control of andarine must be performed in urine on its two major metabolites which were much longer detected than the native molecule.

Novel Aspect

Combination of in-vitro and in-vivo approaches to study by LC-HRMS2 horse metabolism and kinetic of elimination in plasma and urine of andarine in the context of horse doping.

References

1. Hansson A., Knych H., Stanley S., Thevis M., Bondesson U., Hedeland M., Drug Test Anal., 7, 673-83 (2015)

WPS-S06 - MS imaging: applications

Chairs: Martina Marchetti-Deschmann, Fulvio Magni

WP-122 / MICROBIAL MALDI-(2)-MS IMAGING OF STEAM-INACTIVATED CO-CULTURED BACTERIAL COLONIES AT HIGH MASS RESOLUTION

Eike Ulrich Brockmann⁽¹⁾ - ***Daniel Steil***⁽²⁾ - ***Jens Soltwisch***⁽¹⁾ - ***Klaus Dreisewerd***⁽¹⁾

University of Muenster, Institute for Hygiene; Interdisciplinary Center for Clinical Research (IZKF), Münster
(1) - University of Muenster, Institute for Hygiene, Münster (2)

Keywords: microbial MALDI-MSI, MALDI-2, high resolution, pathogenic bacteria, sample preparation

Introduction

Standard protocols to visualize the metabolic exchange in microbial communities by MALDI-MSI are based on analyzing colonies directly from agar [1]. However, this comes with notable drawbacks such as a lacking inactivation step and a limited reproducibility [2]. Here, we introduce a robust membrane-based workflow that enables steam-inactivation and highly-resolved imaging. Examples of competing strains including probiotics (*e.g.* *E. coli* Nissle) are shown.

Methods

E. coli, *S. aureus*, *P. aeruginosa* and *B. subtilis* were investigated as prominent gram-negative and gram-positive samples. Bacterial colonies were co-cultured on LB agar, using mixed ester cellulose membrane filters. Steam inactivation was performed in a home-made chamber. Synapt G2-S QTOF and QExactive plus orbitrap mass spectrometers, equipped with MALDI-2 ion sources (from Waters [3] and Spectrograph [4], respectively), were used for data acquisition.

Results

Inactivation of pathogenic bacteria and metabolic quenching is an indispensable prerequisite for a wider use of microbial MALDI-MSI. With our protocol, this requirement was – for non-sporulating strains – achieved by steam inactivating of membrane-grown colonies for 5 s. Following heat-induced dehydration, the 3-dimensionality of the cultures was reduced to a 2D space. Following matrix application, sample flatness fostered reproducible MALDI analyses and use of laser-induced post-ionization (MALDI-2), boosting the signals of numerous metabolites in positive and negative ion mode. Differentiation and identification of chemical identities was supported by the high resolving power (280k @ m/z 200) and tandem MS capability of the QExactive plus.

Making use of these advantages, we recorded metabolic profiles from several single and co-cultured microbial communities at high lateral and high mass resolution. Examples for the microbial interaction on a molecular level of probiotic *E. coli* Nissle and pathogenic *E. coli* will be demonstrated.

Conclusions

With the emerging interest in microbial MALDI-MSI, a reliable protocol for handling of pathogenic samples is still challenging. Our new membrane-based culturing protocol solves this issue for non-sporulating gram-negative and gram-positive bacteria. The use of MALDI-2 in combination with a high-resolving orbitrap increased chemical coverage, a critical step for a better understanding of the complex chemical communication in bacterial communities.

Novel Aspect

Advanced protocol for microbial MALDI- and MALDI-2 imaging, visualizing the metabolic exchange in (co-)cultured colonies with high sensitivity and mass accuracy.

References

1. Yang Y. J., Phelan V. V., Simkovsky R., Watrous J. D., Trial R. M., Fleming T. C., Wenter R., Moore B. S., Golden S. S., Pogliano K., Dorrestein P. C., *J. Bacteriol*, 194, 6023-6028 (2012)
2. Dunham S. J. B., Ellis J. F., Li B., Sweedler J. V., *Acc. Chem. Res.*, 50, 96-104 (2017)
3. Soltwisch J., Ketting H., Vens-Cappell S., Wiegelmann M., MÜthing J., Dreisewerd K., *Science*, 348, 211-215 (2015)
4. Ellis S. R., Soltwisch J., Paine M. R. L., Dreisewerd K., Heeren R. M. A., *ChemCommun*, 53, 7246-7249 (2017)

WP-123 / BIOIMAGING OF RAT KIDNEYS AFTER SINGLE ADMINISTRATION OF IODINATED CONTRAST AGENTS (ICAS)

***Corinna Brungs*⁽¹⁾ - *Maria Viehoff*⁽¹⁾ - *Robin Schmid*⁽¹⁾ - *Michael Sperling*⁽²⁾ - *Uwe Karst*⁽¹⁾**

Institute of Inorganic and Analytical Chemistry, University of Münster, Münster (1) - European Virtual Institute for Speciation Analysis, University of Münster, Münster (2)

Keywords: iodinated contrast agents, bioimaging, LA-ICP-MS, μ XRF, MALDI-MS

Introduction

Iodinated contrast agents (ICAs), *e.g.*, iopromide and iodixanol, are excreted via the kidneys. The fast and complete excretion of the unmetabolized CA is crucial with respect to toxicity issues. Rats were treated with iopromide or iodixanol with a dosage of 4 g iodine/kg bodyweight. They were sacrificed after 24 or 48 hours, kidneys were embedded and cut into 10 μ m cryosections. To investigate the accumulation and clearance of the ICAs, different bioimaging methods were used.

Methods

Micro X-ray fluorescence (μ XRF) and laser ablation-inductively coupled plasma-mass spectrometry (LA-ICP-MS) were used to investigate the elemental distribution, especially of iodine, in rat kidneys. Matrix-assisted laser desorption/ionization-MS (MALDI-MS) generated information on the distribution of different small molecules in the range of m/z 100-1600 with MS/MS identification.

Results

Significant differences in the accumulation and distribution of iopromide and iodixanol were detected with the three imaging methods when comparing kidney sections of rats, which have been treated with either ICA. In kidney sections treated with iodixanol and analyzed with LA-ICP-MS, the iodine was concentrated in the cortex region. The signal intensities decreased after 48 h. For iopromide, the iodine was accumulated in the medulla region. It was possible to quantify the different concentrations in kidneys spatially resolved. With μ XRF the same distribution and the decrease in intensity after 48 h were investigated.

Iodixanol was detected and identified with MS/MS via MALDI-MS. Additional information about small molecules were possible to analyze and to visualize different regions in tissue sections.

Conclusions

It is possible to investigate the difference in clearance and accumulation in rat kidneys with the different bioimaging methods. Due to the lower iodine concentration in medulla or pelvis, iopromide has a faster clearance. Iodixanol causes a higher urine viscosity, accordingly the iodine amount is higher and it accumulates in the cortex region. Additionally, MALDI-MS can be used to determine the phospholipid distribution in different regions of kidneys.

Novel Aspect

Three different bioimaging techniques are used to investigate the distribution of iopromide and iodixanol in rat kidneys.

WP-124 / REVERSING THE POLARITY – TARGETED MASS-SPECTROMETRY IMAGING ANALYSIS OF FATTY ACIDS IN POSITIVE DETECTION MODE

Andreas Dannhorn⁽¹⁾ - **Emine Kazanc**⁽¹⁾ - **Paolo Inglese**⁽¹⁾ - **James McKenzie**⁽¹⁾ - **Nicole Strittmatter**⁽²⁾ - **Richard J. A. Goodwin**⁽²⁾ - **Zoltán Takáts**⁽¹⁾

Imperial College London, Department of Surgery and Cancer, London (1) - AstraZeneca, Pathology, Drug Safety and Metabolism, IMED Biotech Unit, Cambridge (2)

Keywords: Fatty acids, Charge-remote-fragmentation, metal adducts

Introduction

Qualitative and quantitative analysis of fatty acids can be difficult based on lack of fragmentation shown by

$[M-H]^-$ species. Ionization as metal adducts has proven to be a powerful tool for structure elucidation. Induction of charge-remote fragmentation results in generation of meaningful fragments, enabling identification and quantification of structural isomers.

Methods

To study desorption and fragmentation behavior desorption electrospray ionization (DESI) analysis was performed on mixtures of fatty acids and metals spotted onto PTFE targets. Tandem MS experiments were performed to identify the underlying fragmentation mechanistic. Suited metals were evaluated for direct ionization of fatty acids adducts directly from tissue sections and their fragment ions were used to create a targeted imaging approach.

Results

Lithium adducts showed mixed fragmentation pattern along the fatty acid chain including formation of radical ions while silver adducts showed a preferred cleavage in allyl position to C=C bonds through formation of cyclic intermediates.

The identified fragmentation pattern were subsequently used to identify fatty acids present in tissue sections. Washing steps to remove salts from the tissue sections prior spray-deposition of salt solutions improved the sensitivity of the method. Fatty acids present in tissue sections were identified by DESI-MS/MS experiments. Unique fragments enabling differentiation of structural isomers, *e.g.* arachidonic acid and eicosatetraenoic acid, were subsequently used to create multiple-reaction-monitoring (MRM) transitions for a targeted imaging approach performed on a triple-quadrupole mass spectrometer. The methodology was then used to map the spatial distribution of structural fatty acid isomers across various rodent tissues and human tumor sections.

Conclusions

Ionization of fatty acids as metal adducts has been demonstrated to be a powerful tool for structural identification of fatty acids. Incorporated into a targeted MSI approach it enables to study distribution of fatty acids and potentially of closely related molecules, *e.g.* eicosanoids and prostaglandins in tissues. This is particularly interesting for the study of inflamed tissues and tumors studying the underlying disease mechanisms involving free fatty acids.

Novel Aspect

In this study we present the use of charge-remote fragmentation for identification and targeted analysis of fatty acids in MSI experiments.

WP-126 / MASS SPECTROMETRY IMAGING IN SKIN WITH AP-MALDI AND TOFSIMS

Dana El Assad - Gilles Frache

Luxembourg Institute of Science and Technology (LIST), Materials Research and Technology, Belvaux

Keywords: Imaging mass spectrometry, skin

Introduction

Latest developments in the field of MALDI imaging has led to a significant impact in the pharmaceutical and cosmetics field. MALDI imaging of biological tissue requires high selectivity which is provided by tandem MS and/or accurate mass measurements. Highly selective MALDI imaging is now enabled on high resolution Orbitrap instruments by means of an Atmospheric Pressure MALDI source.

Methods

An AP-MALDI UHR ion source (Masstech Inc.) coupled to a LTQ/Orbitrap Elite high resolution mass spectrometer (Thermo Scientific) was used for the development of targeted and untargeted imaging Mass Spectrometry experiments. Cryosections of biological tissue were coated with matrices using a HTX-TM

sprayer. Time-Of-Flight Secondary Ion Mass Spectrometry (TOFSIMS) was used as a complementary MS imaging modality.

Results

Several applications of the AP-MALDI Orbitrap setup are presented based on samples from various types and sizes for the detection of targeted and untargeted molecules. The instrument allows for the localization of endogenous or exogenous molecules in various types of tissue.

MALDI images are correlated with other imaging techniques in a multimodal imaging workflow including *e.g.* optical microscopy and TOFSIMS.

Conclusions

We show the characteristics and applications of a Masstech AP-MALDI UHR ion source as a part of a multimodal imaging workflow. With a laser spot size down to 10 micron lateral resolution, when associated to a LTQ/Orbitrap mass spectrometer.

Novel Aspect

Multi-modal MS imaging workflow including MALDI HR-MS or -MS_n imaging capabilities down to 10 micron lateral resolution.

WP-128 / DESI MASS SPECTROMETRY IMAGING OF THE DISTRIBUTION OF A FLUORESCENCE MARKER IN MALIGNANT GLIOBLASTOMA BRAIN TUMOURS

Matthew Gentry⁽¹⁾ - ***Emrys Jones***⁽²⁾ - ***David Coope***⁽³⁾ - ***Adam McMahon***⁽¹⁾

University of Manchester, Division of Informatics, Imaging and Data Sciences, Manchester (1) - Waters Corporation, R&D Applications, Wilmslow (2) - Salford Royal NHS Foundation Trust, Neurosurgery, Salford (3)

Keywords: Glioblastoma, DESI-MS, Fluorescence, Imaging, MSI

Introduction

The use of fluorescence guided surgery (FGS), using 5-aminolevulinic acid (5-ALA), is important to establish safe maximal resection and better progression-free survival for glioblastoma (GBM) patients. 5-ALA acts as an orally administered fluorescent pro-agent of protoporphyrin IX (PpIX) (1). Accumulating in glioma cells, 5-ALA is metabolized to PpIX (2), affording intraoperative fluorescence visualization of the tumour, (λ_{ex} 405 nm - blue, λ_{em} 705 nm - red).

Methods

The distribution of PpIX was investigated using Desorption Electrospray Ionization MS Imaging (DESI-MSI) in preclinical and human GBM tissue samples. The fluorescence signature of PpIX was obtained using a fluorescence microscope (Olympus BX51 with QImaging Retiga 6000 camera) with excitation at 405 nm. Mass spectrometry images of the PpIX; parent ion at m/z 563.26 and CID MS/MS fragments, were obtained on a Waters Xevo G2-XS Q-ToF in positive ion mode.

Results

The distribution and detection levels of PpIX in human glioblastoma tissue is reported, with the goal of understanding where 5-ALA, and consequently PpIX, is localised relative to the tumour margins. DESI-MSI analysis of various PpIX dosed tissue homogenates, reveal detection of PpIX at the same concentrations required to observe fluorescence during surgery. Further, we describe the distribution of PpIX in samples from patients who have undergone tumour resection following the administration of 5-ALA at 20 mg/kg, 2 hours prior to time of surgery. DESI-MS images and fluorescence images of PpIX are co-registered alongside

Hematoxylin and Eosin (H&E) staining sections to determine the effectiveness of fluorescence to define tumor margins, with mass spectrometry confirming the molecular source of the fluorescent signal.

Conclusions

The preclinical data confirm that PpIX, the fluorescent porphyrin that results from 5-ALA administration, is detectable by DESI-MSI. By co-registration of fluorescence and mass spectrometry images, we demonstrate that PpIX is detectable to the same concentrations required for fluorescent imaging. This work also measures the relative co-localization of molecular PpIX with the fluorescence signal in the tumours and its ability to define tumour margins.

Novel Aspect

For the first time, the distribution of a fluorescence guided surgical molecule has been mapped in malignant glioblastoma tumour samples using MSI approaches.

References

1. Stummer W., et al., *The Lancet Oncology*, 7, 392-401 (2006).
2. Stummer W., et al., *Journal of photochemistry and photobiology B: Biology*, 45, 160-169 (1998).

WP-129 / MALEIC ANHYDRIDE PROTON SPONGE AND ITS DERIVATIVES TOWARD MALDI-MSI BASED BRAIN TUMOR DIAGNOSIS.

Marco Giampà⁽¹⁾ - ***Annika Zurowietz***⁽¹⁾ - ***Tobias Lammert***⁽¹⁾ - ***Jens Fuchser***⁽²⁾ - ***Volkmar H. Hans***⁽³⁾ - ***Oliver Gembruch***⁽⁴⁾ - ***Udo Kellner***⁽⁵⁾ - ***Tim W. Nattkemper***⁽⁶⁾ - ***Hanna Bednarz***⁽¹⁾ - ***Karsten Niehaus***⁽¹⁾

Bielefeld University, Proteome and metabolome research, Bielefeld (1) - Bruker Daltonics GmbH, BU Pharma, Bielefeld (2) - Universität Bochum, Institut für Pathologie Ruhr, Bielefeld (3) - Universitätsklinikum Essen (AöR), Klinik für Neurochirurgie, Bielefeld (4) - Johannes Wesling Klinikum Minden, Institut für Pathologie, Minden (5) - Bielefeld University, Biodata Mining Group, Bielefeld (6)

Keywords: MALDI-MSI, Cancer diagnosis, Method development, Organic synthesis, Metabolomics

Introduction

The diagnosis of malignant glioma, one of most aggressive brain tumor, is still a challenge in the biomedical and surgical imaging procedures. An emerging cancer hallmark is the reprogramming metabolism where the energy metabolism is adjusted to enhance the energy supply and therefore the biosynthetic pathways for the uncontrolled cell growth. The main metabolic phenotypes in a tumor cell are an overproduction of lactate, so called Warburg effect, and several mutations of mitochondrial enzymes. Due to this, an alteration of cellular metabolism occurs causing the cancer progression [1]. Imaging of cancer metabolism is an optimal diagnostic tool for an accurate glioma classification and localization of tumor in surgical biopsies.

In this regard, Matrix Assisted Laser Desorption Ionization Mass Spectrometry Imaging (MALDI-MSI), a multi-component imaging technique able to detect and localize onto the tissue many compounds at the same run, could perform a visualization of the cancer metabolism directly on the glioma tissue. Unfortunately the detection of metabolites is still challenge because of low concentrations, matrix interference and fragmentations phenomena. It is known that the chemistry of matrix serves fundamental functions in desorption and ionization. For this reason the search for new MALDI matrices is required and then improve the limit detection of metabolites without any matrix background. For the first time, the 1,8-bis(dimethylamino)naphthalene (DMAN) or Proton Sponge was used in MALDI-MS for the detection of small weak acidic compounds with any matrix interference [2]. Unfortunately it was observed that DMAN is vacuum stable and then not suitable for the long-time MALDI-MSI procedures. Indeed in order to increase the vacuum stability the DMAN was structural modified generating two derivatives 1,8-Di(piperidinyl)-

naphthalene (DPN) [3] and 3-(4,5-bis(dimethylamino)naphthalen-1-yl)furan-2,5-dione (4-maleicoanhydrido proton sponge or MAPS) [4],[5].

Regarding the para-substituted derivative MAPS, it was shown that this novel MALDI-MSI matrix is able to detect and localize compounds as small as chloride directly onto glioma tissues [6].

This work wants to show the potentiality of the MAPS as MALDI-MSI matrix in the discovery of novel biomarker candidates in the glioma diagnosis and as synthetic precursor of a new proton sponge based matrices generation.

Methods

Fresh-frozen glioma sections of five patients were covered with a matrix suitable for the detection of small molecules, Maleic Anhydride Proton Sponge (MAPS), using Image Prep™ as sprayer. Afterwards the glioma sections were analyzed by MALDI-ToF (ultrafleXtrem; Bruker) using a spatial resolution of 70 μm. The ionic spatial signatures on the tissue were generated using FlexImaging (Bruker) software. The chemical identification of ions was performed using High resolution MALDI-FT-ICR (solarix; Bruker) Imaging. Furthermore, in order to investigate the tissue microheterogeneity, clustering analysis of MALDI-MSI data was performed using Cardinal® as statistical software [6]. The cluster analysis was compared with histological characterization of H&E stained glioma sections, using QuPath as software.

Two proton sponge derivatives, 4-(*N*-phenyl)maleimido proton sponge (PMPS) and 4-(*N*-methyl)maleimido proton sponge (MMPS), were synthesized and sprayed onto mouse brain sections using an airbrush. MALDI-MSI experiments were performed in reflectron negative mode using 100 μm of spatial resolution and compared with MAPS.

Results

Several metabolic pathways of glioma metabolism, *e.g.* aerobic glycolysis, tricarboxylic acid cycle, and aminoacids metabolism, were visualized by MALDI-MSI procedures.

In the aerobic glycolysis, ionic maps of intermediates such as glucose, glucose 6-phosphate, and lactate were observed, confirming the Warburg effect as general biomarker of cancer tissue.

MALDI-MSI of citric acid cycle metabolites, succinate and fumarate, were detected and different spatial signatures were observed in five patients thus suggesting mitochondrial enzymatic aberrations. It is very interesting to note that the spatial distribution on the tumor region of amino acids such as glutamine, which has a crucial role in the cancer progression. Moreover based on metabolite's distribution and through an unsupervised cluster analysis, specific clusters were found in the five patients showing a correlation with tumor type and grading. Furthermore the intra-tumor differentiation, shown by the clustering analysis, was confirmed by PCA.

For the first time the 4-(*N*-phenyl) proton sponge (PMPS) and 4-(*N*-methyl) proton sponge (MMPS) were synthesized, characterized and applied as matrix in the MALDI-MSI of mouse brain tissue section. Small metabolites were detected and imaged on the tissue. Indeed several tissue-specific ionic maps were obtained showing a correlation with the histological features.

Conclusions

MALDI-MSI of glioma metabolism showed the detection of patient survival correlated biomarkers by phenotypically spatial signatures *e.g.* IDH1 or IDH2 mutations. Moreover it was shown the heterogeneity in metabolite distribution, confirming the tissue microheterogeneity that is one of main reasons of drug resistance. Moreover a new generation of proton sponge were synthesized and applied in MALDI-MSI showing very good performance as matrices. In conclusion the para-substituted proton sponges are valid candidate as matrix for tissue profiling by using MALDI-MSI procedures.

Novel Aspect

MAPS and its derivatives could contribute to the cancer diagnosis by MALDI-MSI procedures. Moreover the generation of new proton sponge derivatives could permits to gain more insight into the MALDI of metabolites. In conclusion this study is an emerging contribution for the clinical tumor diagnosis and also an outlook for a more detailed investigation of the ion formation mechanisms in MALDI that occur in complex systems such as biological tissue.

References

1. Masui, K., Cavenee, W. K., Mischel, P. S., Brain tumor pathology, 33, 161-168 (2016).
2. Shroff, R., Svatoš, A., Analytical chemistry, 81, 7954-7959 (2009).
3. Weißflog, J., Svatoš, A. RSC Advances, 6, 75073-75081 (2016).
4. Swor, C. D., Zakharov, L. N., Tyler, D. R. The Journal of organic chemistry, 75(20), 6977-6979 (2010).
5. Giampà, M., Lissel, M. B., Patschkowski, T., Fuchser, J., Hans, V. H., Gembruch, O., Bednarz H. Niehaus, K., Chemical Communications, 52, 9801-9804 (2016).
6. Bemis, K. D., Harry, A., Eberlin, L. S., Ferreira, C., van de Ven, S. M., Mallick, P., Stolovitz M., Vitek, O., Bioinformatics, 31, 2418-2420 (2015).

WP-130 / THE EFFECT OF PROTEIN BINDING OF DRUG SUBSTANCES IN DESORPTION ELECTROSPRAY IONIZATION MASS SPECTROMETRY IMAGING

Anne Mette Handler⁽¹⁾ – *Kim Troensegaard Nielsen*⁽²⁾ – *Christian Janfelt*⁽²⁾

LEO Pharma, University of Copenhagen, Copenhagen (1) - University of Copenhagen, Department of Pharmacy, Copenhagen (2)

Keywords: Protein-drug interactions, Mass Spectrometry Imaging, Desorption Electrospray Ionization

Introduction

Mass Spectrometry Imaging (MSI) is increasingly used as a tool in pharmaceutical R&D for analysing the distributions of drugs in tissue [1]. However, some drug substances are detected with greatly compromised sensitivity compared to traditional bioanalysis, due to strong binding to blood proteins such as albumin [2,3]. The study characterizes the protein binding of drugs in MSI and different approaches to improve the MSI detection of protein bound drug substances.

Methods

Cryo-sections were made of liver homogenates spiked with different drug substances, representing different degrees of protein binding. For each cryo-section, one half was treated with different denaturants in order to break the protein binding. The sections were analyzed by Desorption Electrospray Ionization (DESI)-MSI, and the treated and non-treated tissue was compared.

Results

DESI-MSI was applied to characterize protein binding, and significantly reduced signals were observed in tissue homogenates spiked with a strongly protein bound drug substance like Roflumilast, compared to cryo-sections of blocks of gelatin spiked to the same concentration. This effect was less pronounced for drug substances like codeine known to have a lower degree of protein binding. Different denaturants are being tested on similar cryo-section in attempts to alleviate the protein binding and improve the DESI-MSI sensitivity (work still in progress). Some of the denaturants appear to be more efficient than others in breaking drug-protein bonds.

Conclusions

The use of denaturants provide an increase in the free fraction of drug which can be detected by DESI-MSI. Although the results were obtained by DESI-MSI, it is likely that similar findings would be made with other MSI methods such as MALDI-MSI. The use of agents that denature proteins thereby breaking the drug-protein binding is worth considering as a sample preparation step in MSI in imaging of drug substances which are known to have strong protein binding.

Novel Aspect

Investigation of protein-drug interactions in Mass Spectrometry Imaging and methods for disruption of the

protein binding.

References

1. Swales J., Hamm G., Clench M., Goodwin R., International Journal of Mass Spectrometry, In Press, Available online: March 2018
2. Chumbley, C. W., Reyzer M. L., Allen J. L., Marriner G. A., Via L. E., Barry C. E., Caprioli R. M., Analytical Chemistry, 88, 2392-8 (2016)
3. Bonnel D., Legouffe R., Eriksson A. H., Mortensen R. W., Pamelard F., Stauber J., Nielsen K. T., Analytical and Bioanalytical Chemistry, 410, 2815-2828 (2018).

WP-131 / MS FOR MS; A MASS SPECTROMETRY-BASED STUDY INTO MULTIPLE SCLEROSIS

Bram Heijs⁽¹⁾ - ***Floris-Jan de Graaff***⁽¹⁾ - ***Evelyne Steenvoorden***⁽¹⁾ - ***Martin Giera***⁽¹⁾ - ***Gijs Kooij***⁽²⁾

Center for Proteomics & Metabolomics, Leiden University Medical Center, Leiden (NL) (1) - 2Department of Molecular Cell Biology and Immunology, MS Center Amsterdam, Neuroscience Campus Amsterdam, VU University Medical Center, Amsterdam (NL) (2)

Keywords: Mass spectrometry imaging; quantitative LC-MS/MS; Lipidyzer; Multiple Sclerosis; Lipids

Introduction

Multiple Sclerosis (MS) is a chronic inflammatory disease of the central nervous system (CNS) in which infiltrating immune cells form MS lesions, resulting in extensive demyelination, axonal damage and neurodegeneration. Lipids are important mediators in immune response and a main constituent of the CNS, but their phenotype and contribution to MS pathogenesis remain elusive. Here we study lipids by applying multiple mass spectrometry based approaches to human MS.

Methods

Human post-mortem brain tissue from 3MS patients and two healthy controls were collected and sectioned for both high resolution MALDI-TOF-MSI and MALDI-FTICR-MSI, as well as for quantitative LC-DMS-MS/MS. To confirm the location of MS lesions, we performed immunohistochemical staining (IHC), using anti-myelin (PLP) and immune cell markers (LN3). The MSI data was interrogated in a histology-directed manner, and by using multivariate statistical methods.

Results

Interrogation of the quantitative LC-MS/MS data, acquired using the SCIEX Lipidyzer platform after homogenizing tissue sections, resulted in the finding that several ceramides (dLCER(16:0), LCER(18:1), LCER(18:0)) were considerably higher in concentration in the MS patient tissues, compared to the healthy control tissues. Based on mass matching, *m/z* features corresponding to different ceramides, both lactosyl-, and hexosylceramides (HCER), were detected specifically in the MS lesions in the MALDI-MSI analyses. The specific detection of overexpressed LCER lipids in MS confirms the findings from a previous study on an experimental autoimmune encephalomyelitis (EAE) mouse model of MS, where astrocyte-derived LCER lipids were shown to contribute to the EAE disease phenotype [1]. Besides the LCER and HCER lipids, additional lipid features were found to be specific to the MS lesions, and using spatial segmentation of the tissues the MS lesions were clearly discerned from surrounding tissues.

Conclusions

In this proof-of-principle study on the investigation of the MS lipidome by using several MS techniques, we were able to confirm, for the first time, the presence and increased abundance of lactosylceramides inside the lesion of human MS patients. This study provides a solid basis into the further investigation of the role and location as well as function of these lipids in the context of MS pathogenesis.

Novel Aspect

The combination of off-tissue quantitative LC-DMS-MS/MS and in-situ high-resolution MALDI-MSI for the study of lipids in human MS samples.

References

1. Rostami A., Ciric B., Nature Medicine, 20, 1092-1093 (2014)

WP-132 / GOLD-COATED BLACK SILICON: AN EFFICIENT SUBSTRATE FOR LASER DESORPTION IONIZATION MASS SPECTROMETRY

Stefania-Alexandra Iakab - Pere Rafols - Oscar Yanes - Xavier Correig

Universitat Rovira i Virgili, Department of Electrical Electronic Engineering and Automation, Tarragona

Keywords: black silicon, Au, LDI-MSI, tissue imprinting

Introduction

Silicon based substrates such as DIOS and NIMS proved their value in MS studies by successful detection of small molecules with reduced background noise and fragmentation, despite their unsafe fabrication method [1,2]. Au nanoparticles were used in laser desorption ionization MS analyses with increased resolution and sensitivity [3,4]. We developed Au coated black silicon substrates (BSi) with hydrophobic/hydrophilic selective regions for MSI applications.

Methods

The BSi substrates were created by dry etching with SF₆ and O₂ plasma mixture [5]. Au layer was deposited by sputtering and the hydrophobic or hydrophilic surfaces were created with CHF₃ or O₂ plasma, respectively. Substrates were analysed by SEM, AFM, UV-VIS spectroscopy and mean contact angle. Animal tissue samples were imprinted on the substrates by touching. The MSI acquisition was executed with a MALDI TOF/TOF UltrafleXtreme spectrometer from Bruker Daltonics.

Results

Physicochemical characterization of the AuBSi substrates showed a needle-like array structure (<900 nm height) that absorbs ~95% of visible light, specifically 98% at 355 nm (the wavelength of the MALDI Nd:Yag laser). The hydrophobic regions of the substrates have a mean contact angle (CA) of ~165°, because of the CH_x terminal groups, while the hydrophilic regions have a mean CA of ~45°, due to the HO groups. The AuBSi substrates revealed their efficacy by promoting the ionization/desorption processes of standard compounds as well as metabolites in tissue samples. The detected Au cluster peaks allow an accurate mass calibration process throughout the full mass range. Many imprinted endogenous molecules were detected in the MS images. Therefore, the selectivity of the substrate was demonstrated by the specific adhesion of molecules: hydrophobic or hydrophilic compounds adhered to hydrophobic or hydrophilic regions, respectively. This may simplify untargeted analysis (molecule assignment process) of polar analytes from a broad mass range.

Conclusions

Selective imaging of adhered compounds from imprinted tissues is possible with our novel selective AuBSi substrates. Compounds adhered to the AuBSi substrate based on polar interactions between tissue and surface: hydrophobic analytes were found in the hydrophobic region and hydrophilic analytes were found in the hydrophilic region. Thus, knowing the polarity of the detected molecules simplifies the putative identification process.

Novel Aspect

For the first time selective AuBSi substrates were entirely fabricated by dry etching and sputtering. No

sample preparation was needed for selective imprinting of molecules.

References

1. Wei, J., Buriak, J. M., Siuzdak, G., Nature 399, 243-246 (1999).
2. Northen, T. R. et al., Nature 449, 1033-1036 (2007).
3. McDonnell, L. A., Heeren, R. M. A., Mass Spectrom. Rev. 26, 606-643 (2007).
4. Chughtai, K., Heeren, R. M. A., Chem. Rev. 110, 3237-3277 (2011).
5. Gao, J., De Raad, M., Bowen, B. P., Zuckermann, R. N., Northen, T. R., Anal. Chem, 88, 1625-1630, (2016).

WP-133 / APPLICATION OF ISOTOPE-LABELING IN WHOLE-BODY DESI-MSI STUDIES FOR INTRA-ANIMAL COMPARISON OF DRUG ADMINISTRATION ROUTES

Loïs Morineau⁽¹⁾ - *Sophie Chakroun Jacobsen*⁽²⁾ - *Karen Kleberg*⁽³⁾ - *Harald S. Hansen*⁽³⁾ – *Christian Janfelt*⁽²⁾

Université de Nantes, UFR Sciences et Techniques, Nantes (1) - University of Copenhagen, Department of Pharmacy, Copenhagen (2) - University of Copenhagen, Department of Drug Design and Pharmacology, Copenhagen (3)

Keywords: Drug distribution analysis, DESI-MSI, whole-body imaging, stable isotope labeling, drug administration

Introduction

The distribution of an administered drug in the organism is highly dependent on the administration route as well as the drug formulation. Therefore, assessing the distribution and metabolism of drug substances or drug candidates is of great importance in drug development, and mass spectrometry imaging has demonstrated many advantages, in particular due to its unique selectivity which enables specific imaging of several exogenous and endogenous compounds simultaneously.

Methods

The antidepressant drug amitriptyline was intravenously (IV) dosed to a mouse and immediately after, a similar amount of a deuterium labeled version of the drug was intraperitoneally (IP) dosed to the same animal. Whole-body cryo-sections were made at t=5, 15, 30 and 60 min, according to a previously published protocol [1], and the two drug substances were simultaneously imaged with high sensitivity by DESI-MS/MS using the Displaced Dual-mode imaging approach [2].

Results

The approach provides a unique way to compare administration routes in the very same animal, thereby eliminating uncertainties caused by animal-to-animal variation. Already after 5 minutes, the IV dosed drug was detected throughout the animal, while the IP dosed drug (deuterium labeled) was primarily found in the abdominal cavity. After 15 minutes, the IP injected drug appeared highly similar to the 5 minutes IV injection, while the IV injected had started accumulation in the stomach region. At later survival times, the differences of the two administration routes became less pronounced.

The deuterium label was partly preserved in metabolism to the primary metabolite nortriptyline, and the approach thus also enabled assessment on how the drug metabolism depends on the administration route. Interestingly, even at early survival times, the two administration routes provided highly similar metabolite distributions.

Conclusions

The presented method does not only enable IV and IP comparison, but could also be used for comparison of other administration routes such as subcutaneous or oral administration, and in the case of the latter,

two different oral formulations could be tested in the same animal. While autoradiography has higher sensitivity, this study demonstrates how much more information MSI can provide with its unique selectivity, even within different isotopologues of the same drug substance.

Novel Aspect

Highly selective MS imaging of different isotopologues of a drug and its metabolites provides unprecedented Intra-animal comparison of administration routes.

References

1. Okutan, S. ; Hansen, H. S. ; Janfelt, C. *Proteomics* 2016, 16, 1634–1641.
2. Janfelt, C. ; Wellner, N. ; Hansen, H. S. ; Hansen, S. H. *J. Mass Spectrom.* 2013, 48, 361-366.

WP-134 / SINGLE CELL LIVESCAN MASS SPECTROMETRY PROFILING AND IMAGING

Filip Kaftan⁽¹⁾ - **Marine Vallet**⁽²⁾ - **Tim Baumeister**⁽²⁾ - **Georg Pohnert**⁽²⁾ - **Aleš Svatoš**⁽¹⁾

Max Planck Institute for Chemical Ecology, Research Group Mass Spectrometry/Proteomics, Jena (1) - Max Planck Institute for Chemical Ecology, Research Group Plankton Community Interaction, Jena (2)

Keywords: LiveScan, single cell MS/MSI, diatoms, molecular imaging, surface associated metabolites

Introduction

In recent years, there has been a growing interest in the investigation of single cells and unicellular organisms. Single-cell measurement science has numerous applications across the biomedical and life science fields. Especially development of ambient MS techniques along with high performance MS systems provides powerful tools [1] which push the boundaries of system sensitivity so that only a minimal amount of cell material is needed to perform an analysis.

Methods

Unicellular marine diatoms (*Coscinodiscus* sp.), freshwater algae (*Haematococcus* sp.) and various plant cells were analyzed in LDI mode via AP-SMALDI10-Orbitrap-MS/MSI Q-ExactivePlus MS system (TransMIT/Thermo Fisher Scientific, Germany). Cells were picked up from their media and transferred on a wet membrane filter or placed in a thin hydrogel film. All substrates were moistened/rehydrated with a suitable medium in order to ensure the cell viability.

Results

The nanostructure of the *Coscinodiscus* sp. silica frustule acts as UV laser energy absorber and provides an intense ion signal in LDI mode [2] without MALDI matrix. Especially pigments characteristic to algae such as chlorophyll a (m/z 892), fucoxanthin (m/z 658) or beta-carotene (m/z 536) were present in the spectra together with their common fragments. The identity of these signals was confirmed with standards. Furthermore amino acids, fatty acids and different groups of lipids were found in the MS spectra. In the same way an AP-SMALDI-FT MS analysis of *Haematococcus pluvialis* cells showed strong signals of its typical red pigment, the antioxidant astaxanthin (m/z 596). Moreover, PLS-D analysis of the MS spectra showed a clear separation between groups of nutrient-depleted cultures of *Haematococcus pluvialis* cells (green phase) and repleted cultures (red phase). Using the very same approach it was also possible to show distinctive patterns of old and young cells of *Coscinodiscus granii*.

Conclusions

Our LiveScan method requires only little sample preparation. Fresh and vital algae cells isolated from their aqueous substrata were placed on a wet membrane filter/hydrogel and analyzed directly. The obtained MS spectra covered different groups of compounds including, amino acids, fatty acids, pigments and potential

signaling molecules in high intensity. Characteristic patterns of the acquired spectra enable the differentiation of cell species.

Novel Aspect

LiveScan MS-MSI represents a simple, fast and high throughput method for analyzing single cells in liquid phase under ambient conditions.

References

1. Khalil S. M., Pretzel J., Becker K., Spengler B., *IJMS* 416, 1 (2017).
2. Jaschinski T., Svatos A., Pohnert G. *Rapid Commun. Mass Spectrom.* 27, 109 (2013).

WP-135 / EVALUATION OF EMBEDDING STRATEGIES FOR FRAGILE AND SMALL SPECIMEN IN MSI ANALYSIS

Emine Kazanc - Andreas Danhorn - Paolo Inglese - James McKenzie - Zoltán Takáts

Department of Surgery and Cancer, Sir Alexander Fleming Building, Imperial College London, London SW7 2AZ, UK.

Keywords: MSI, fragile specimen, tissue embedding

Introduction

Sectioning of fragile and small specimen can be challenge for mass spectroscopy imaging(MSI) analysis. However, embedding multiple small specimen allows simultaneous sectioning. An optimum embedding medium should meet several requirements, *e.g.* no altering the tissues, no negative effects for MSI and be sectionable of multiple specimen. The comparison of polymers with the freezing methods for embedding are critical factors in obtaining the best MSI result.

Methods

5 different polymers (Gelatine, *N*-(2-Hydroxypropyl)-methacrylamide (HPMA), Sodium-Carboxymethyl-cellulose (CMC), (Hydroxypropyl-)methylcellulose (HPMC), Polyvinylpyrrolidone (PVP) were evaluated for their MSI suitability. All embedding media were evaluated for ease of sectioning, adherence to tissues and matrix interference, analyte delocalization and post-DESI histological evaluation.

Results

Preliminary results displayed that Gelatine, HPMA, CMC, HPMC and PVP modified HPMC solution are able to provide easy sectioning whilst CMC with addition of alcohol and PVP were not suitable as they were not fully frozen under usual sectioning conditions and were not further investigated.

Adherence of the embedding medium to the tissue is important as the tissue sections must not fall out or curl up during sectioning. CMC and PVP modified HPMC are well adapted for this purpose while HPMA and HPMC are mostly suited for larger samples but not for small samples like intestines.

The last significant criteria is matrix interference. Gelatine shows chemical background from *e.g.* choline and arginine and tissue damage from hot gelatine solution. CMC showed adduct shift towards Na-adducts and smearing. HPMA, HPMC and PVP modified HPMC showed no matrix interference in MALDI and DESI-MSI.

Conclusions

We found that the PVP modified HPMC is meeting all the criteria of an optimal embedding media, easy sectioning, high adherence to tissues and non-observed matrix interference.

We also found that snap freezing with isopropyl alcohol is a very promising method for preventing analyte delocalisation, block fracturing and small tissues thawing before fully frozen.

Novel Aspect

This study compares different polymers and freezing methods to develop a comprehensive tissue embedding protocol for MSI analysis.

WP-136 / VISUALIZING MOLECULAR DISTRIBUTIONS ACROSS/UNDER THE HUMAN SKIN BY 2D-DESI MS IMAGING OF TAPE-STRIPPED STRATUM CORNEUM

Keishi Kihara - Akira Motoyama

Shiseido, Global Innovation Center, Yokohama

Keywords: DESI, Mass Spectrometry Imaging (MSI), Natural Moisturizing Factor (NMF), Tape-stripping, Stratum Corneum (SC)

Introduction

We have developed a new method combining DESI MS imaging (MSI) [1] and tape-stripping[2] techniques to visualize planer/depth profiles of natural moisturizing factor (NMF). NMF is a mixture of hygroscopic small molecules important to maintain hydration and barrier homeostasis of stratum corneum (SC) [3]. Differences of NMF content/composition are good index of skin conditions, whereas fine molecular distributions could not be obtained by conventional methods.

Methods

SC samples were collected by tape-stripping from the face or the inner forearm of healthy volunteers (one 30s female and one 30s male). Polymer-based adhesive sheets(25 × 70 mm, rectangular, in-house product) were used for SC collection. A high-resolution mass spectrometer equipped with a 2D-DESI source(Synapt-G2-SI, Waters) was used for DESI MSI analysis. Data processing was conducted by HDI-Imaging software® (Waters).

Results

After the careful optimization of DESI experimental conditions, the developed visualizing method allowed to detect important 16 NMF-related molecules simultaneously. As a proof-of-principle experiments, planar distribution of NMF molecules across the half face was visualized. SC samples were collected from one side of the face on 8 adhesive sheets and measured with 500 um spatial resolution. It was found that the distribution/conversion of NMF molecules had large diversity. For example, serine (Ser) and pyrrolidone carboxylic acid (PCA), which are both abundant NMF [3], showed completely different distribution profiles (PCA distributed unevenly opposed to Ser). We also successfully visualized the changes of NMF depth profiles when the skin was treated with 2% sodium lauryl sulfate (SLS). SCs were collected before/after the treatment and measured by DESI-MSI. PCA was significantly eluted by the treatment, while the elution was prevented when a hand cream containing a film-forming ingredient was applied before the SLS treatment.

Conclusions

NMF molecular distributions across the face and depth profile in the arm skin were successfully visualized for the first time. Because of the multi-target simultaneous analysis capability of the developed DESI-MSI method, molecular kinetics such as enzymatic conversion could be assessed. We believe that the new method is a promising tool for skin biological research and will provide useful insights for the development of new skincare products.

Novel Aspect

The first application of DESI MSI to tape-stripped SC to visualize molecular distribution of NMF across the human face and depth profile.

References

1. Takats Z., Wiseman J. M., Gologan B., Cooks R. G., Science, 306, 471 (2004)

- Jacobi U., Weigmann H. J., Ulrich J., Sterry W., Lademann J., *Skin Res. Tech.*, 11, 91 (2005)
- Koyama J., Horii I., Kawasaki K., Nakayama Y., Morikawa Y., Mitsui J., *J. Soc. Cosmet. Chem.*, 35, 183-195 (1984).

WP-137 / DISTRIBUTION OF LIPIDS AND ANTI-TB DRUGS IN MOUSE LUNG TISSUE

Julia Kokesch-Himmelreich⁽¹⁾ - Axel Treu⁽¹⁾ - Kerstin Walter⁽²⁾ - Christoph Hölscher⁽²⁾ - Andreas Römpf⁽¹⁾

Bioanalytical Sciences and Food Analysis, University of Bayreuth, Bayreuth (1) - Infection immunology, Research Center Borstel, Borstel (2)

Keywords: MALDI imaging, anti-TB drugs, high resolution

Introduction

The rise of multi drug-resistant *Mycobacterium tuberculosis* strains creates an ever increasing need for novel anti-TB drugs. The development of anti-TB drugs however is impeded by the shortcomings of current preclinical animal models. In preparation for studies on novel anti-TB drugs, an AP-MALDI imaging method was established in this study to map the distribution of phospholipids and multiple anti-TB drugs in mouse lung tissue.

Methods

MS imaging experiments were carried out on a Q-Exactive™ HF Hybrid-Quadrupole-Orbitrap mass spectrometer (Thermo Fisher Scientific, Bremen), coupled to the AP-SMALDI10 ion source (TransMIT GmbH, Gießen) [1]. Mass spectra with high mass accuracy (<2 ppm) and high mass resolution (R=240,000 @m/z 200) were acquired. Measurements were performed on lung tissue of BALB/c mice treated with anti-TB drugs. Different matrices were applied using a pneumatic sprayer system.

Results

Sample preparation including matrix application and cryosectioning was optimized for several anti-TB drugs. MS imaging measurements of mice lung sections were performed with pixel sizes between 10 and 40 µm. These experiments showed high spatial details in phospholipid images, which corresponded well with histological features. Phospholipid patterns differed substantially between different tissue regions. In addition, the distribution of several anti-TB drugs could be shown in mouse lung tissue.

Conclusions

The investigated anti-TB drugs differ significantly in their physicochemical properties, for example polarity. Therefore, a sectioning protocol was adapted to avoid leaking of the drug from tissue. Different matrices and solvents had to be used for different drugs. Several drugs were detectable in treated mice, but showed different limits of detection. In conclusion, the MALDI imaging workflow must be optimized for each drug.

Novel Aspect

We established a high resolution MALDI MSI method to analyze distributions of lipids and anti-TB drugs in mouse lungs with a spatial resolution close to the cellular level.

References

- Römpf A., Spengler B., *Histochem. Cell Biol.* 139, 759-783 (2013).

WP-138 / CALCIUM AND LIPID IMAGING OF REGENERATED BONY RAYS OF ZEBRAFISH CAUDAL FIN USING TIME-OF-FLIGHT SECONDARY ION MASS SPECTROMETRY

Sun Young Lee - Jae Young Kim - Dae Won Moon

Keywords: ToF-SIMS, MS imaging, zebrafish, bone, small molecule analysis

Introduction

We observed regenerated caudal fins of zebrafish using time-of-flight secondary ion mass spectrometry (ToF-SIMS). The spatial distribution of the small molecules including mineralized ions, lipids and fatty acids in the regenerated area of the fin was found to be different from the spatial distribution in the existing part.

Methods

ToF-SIMS images in this study were obtained using ToF-SIMS 5 (ION-TOF GmbH) equipped with liquid metal ion gun (LMIG) and gas cluster ion beam (GCIB) sources [1,2]. The caudal fins were sputtered with Ar clusters and etched to 20-25 μm from the surface in the depth direction, then ToF-SIMS images of the bony ray and the inter-ray mesenchymal tissue were obtained using rastering of a Bi^{3+} ion.

Results

ToF-SIMS analyzed both positive and negative ion modes and about 1000 specimen-related spectra were obtained from zebrafish caudal fins. Since the ToF-SIMS method detect mineral ions very well, it provided high spatial resolution ion images including Ca^+ , Mg^+ , K^+ , and PO_2^- ions to intuitively identify the precise location of the truncated site and the regeneration shape of the bony rays of the zebrafish caudal fin at 4 days-post-amputation (dpa) and 7 dpa [3]. Whereas the calcium and PO_2^- ions were be found to contain more in the regenerated area than in the existing area of the caudal fin at 4 dpa, there were little difference between the two areas of the caudal fin at 7 dpa. For phosphocholine ion in the caudal fins at 4 dpa and 7 dpa, the regenerated area had a higher signal intensity than the existing area.

Conclusions

Several ion images of regenerated bony rays in the zebrafish caudal fin indicated that the regeneration of bony rays might precede the regeneration of other tissue materials. Consequently, ToF-SIMS imaging is a simple and easy approach to understanding the spatial distribution of small molecules in the regenerated caudal fin of a zebrafish.

Novel Aspect

The zebrafish caudal fin is not only a very appropriate bio-specimen for mass spectrometry due to its thickness of less than 200 micrometers but also a very important organ for tissue regeneration studies.

References

1. Vanbellinghen, Q. P. ; Elie, N. ; Eller, M. J. ; Della-Negra, S. ; Touboul, D. ; Brunelle, A. Rapid Commun. Mass Spectrom, 29, 1187-1195 (2015)
2. Shon, H. K. ; Yoon, S. ; Moon, J. H. ; Lee, T. G. Biointerphases, 11, 02A321 (2016).
3. Pfefferli, C. ; Jaźwińska, A. Regeneration, 2 , 72-83 (2015).

WP-139 / TOF-SIMS IMAGING OF RAT AND DROSOPHILA BRAIN TISSUE SECTIONS PREPARED BY USING THE NEW SAMPLING METHOD

Tae Geol Lee⁽¹⁾ - ***Hyun Kyong Shon***⁽¹⁾ - ***Minh Uyen Thi Le***⁽¹⁾ - ***Jin Gyeong Son***⁽¹⁾ - ***Sohee Yoon***⁽¹⁾ - ***Jeong Hee Moon***⁽²⁾ - ***Chan Young Shin***⁽³⁾ - ***Sung Bae Lee***⁽⁴⁾

KRISS, Center for Nano-Bio Measurement, Daejeon (1) - KRIBB, Disease Target Structure Research Center, Daejeon (2) - Konkuk University, Center of Pharmacology, Seoul (3) - DGIST, Department of Brain and Cognitive Science, Daegu (4)

Keywords: ToF-SIMS imaging, Tape-supporting sampling, Brain, Rat, Drosophila

Introduction

Brain imaging using ToF-SIMS has been reported to produce the distorted biomolecular distributions due to the cholesterol-induced matrix effect, delocalization of the analytes, and/or inadequate tissue adherence during sample preparation. In this talk, we propose new sample preparation and drying methods that can be used to obtain accurate biomolecular images of brains at room temperature, instead of frozen-hydrate analysis using liquid-nitrogen.

Methods

The rat or Drosophila brains prepared by using a double-sided adhesive tape on a pre-cooled ($-20\text{ }^{\circ}\text{C}$) stainless steel plate was freeze-dried in a load-lock chamber of ToF-SIMS for about an hour and moved directly to the main chamber. All measurements (depth profiles, images) were obtained using a TOF-SIMS V instrument. The 25-keV Bi^{3+} (current $\sim 0.1\text{ pA}$, cycle time: $200\text{ }\mu\text{s}$) was used as an analysis ion beam and the 5-keV Ar 1500^{+} (current $\sim 3\text{ nA}$) was used as a sputtering ion beam.

Results

We obtained the depth profiles of rat brain tissue (corpus callosum) as a function of sampling and drying methods. In the white matter of the rat brain prepared by using the thaw-mounting method, cholesterol migrated to the sample surface regardless of room temperature-drying or freeze-drying. However, freeze-drying of the sample prepared by using the tape-supporting method did not cause the migration of cholesterol to the sample surface. In the frozen-hydrate analysis, cholesterol also did not migrate to the sample surface of the rat brain prepared with the thaw-mounting method or tape-supporting method [1]. Regarding the Drosophila brain imaging, the use of conductive double-sided adhesive carbon tape has been shown to enable the successful collection of small and fragile samples to surfaces compatible with ToF-SIMS analysis [2]. In addition, the tape-supporting method can effectively preserve the distribution of biomolecules on the rat and Drosophila brains in the same manner as that of the frozen-hydrated analysis (below -100°C) [1,2].

Conclusions

Accurate and undistorted images of biomolecules in the rat and Drosophila brains can be obtained from brain sections prepared by the tape-supporting and freeze-drying method without performing frozen-hydrate analysis. We expect that this new sample preparation method will enable scientists using ToF-SIMS to easily and accurately detect important biomarkers for degenerative diseases.

Novel Aspect

Our protocol circumvents cholesterol migration to the surface, delocalization of small molecules, and surface water interference.

References

1. Shon H. K., Kim S. H., Yoon S., Shin C. Y., Lee T. G., *Biointerphases*, 13, 03B411 (2018).
2. Le M. U. T., Son J. G., Shon H. K., Park J. H., Lee S. B., Lee T. G., *Biointerphases*, 13, 03B414 (2018).

WP-140 / MALDI MSI AND IMMUNOHISTOCHEMISTRY ANALYSIS OF COLORECTAL CARCINOMA SPHEROIDS

Markéta Machálková⁽¹⁾ - *Jarmila Navrátilová*⁽²⁾ - *Barbora Pavlatovská*⁽²⁾ - *Katarzyna Anna Radaszkiewicz*⁽²⁾ - *Adam Pruška*⁽¹⁾ - *Jan Preisler*⁽¹⁾

Masaryk University, Department of Chemistry, Brno (1) - Masaryk University, Department of Experimental Biology, Brno (2)

Keywords: MALDI, MS imaging, spheroids, confocal microscopy, immunohistochemistry

Introduction

Multimodal imaging is a growing scientific field beneficial for explanation of many biological and chemical processes. Matrix-assisted laser desorption/ionization mass spectrometry imaging (MALDI MSI) together with laser scanning confocal microscopy (LSCM) conducted immunohistochemistry (IHC) allow to study drug distributions and related cellular responses. [1] In our study, 3D cell lines (spheroids) were employed as improved model for drug testing [2].

Methods

Spheroids were treated with the drug perifosine for selected time intervals, transferred to gelatine, frozen, sectioned and placed to conductive MALDI glass slides. Dried sections were covered with sublimated DHB matrix and measured using MALDI MSI. The matrix was washed away by organic solvents and samples were processed using standard IHC protocol with selected antibodies. Later, the analysis on LSCM and its correlation to MALDI MSI was carried out.

Results

The study of perifosine penetration by MALDI MSI revealed that perifosine does not penetrate the spheroid core even after 24 hour induction, it concentrates only in outer 150- μm cell layer. A specific characteristic of this drug is also its diffusion from the tissue to gelatine medium during sample preparation. Viability assays by IHC were conducted for markers of cell proliferation and apoptosis with parallel staining of all cell nuclei. We detected elevated signals of apoptotic marker in spheroid core and proliferative cells on spheroid edges. This corresponds to general knowledge that spheroids from certain size physiologically develop necrotic core, due to by the lack of oxygen and nutrients in deeper cell layers [3]. Unfortunately, none of these spheroid physiologies were changed when the spheroids were treated by perifosine, meaning perifosine does not have any negative effect on cell viability.

Conclusions

MALDI mass spectrometry imaging and immunohistochemistry analysis of colorectal carcinoma spheroids showed no correlation between perifosine and dead cell presence. Apoptotic cells were found in spheroid core, however this was caused only by natural sample physiology. Perifosine was also not capable of penetrating the spheroid deeper than 150 μm , suggesting a poor drug efficacy.

Novel Aspect

The analysis of the same tissue section by both MALDI MSI and immunohistochemistry, enabling correlation of both drug distribution and cell viability.

References

1. Huber, K., et al., Analytical Chemistry, 86(21), 10568-10575 (2014).
2. Fang, Y. and Eglén, R., Slas Discovery, 22(5), 456-472 (2017).
3. Friedrich, J., et al., Nature Protocols, 4(3), 309-324 (2009).

WP-141 / MULTIMODAL IMAGING OF UNDECALCIFIED TISSUE SECTIONS BY MALDI MS AND μXRF

Anastasiya Svirkova⁽¹⁾ - **Anna Turyanskaya**⁽²⁾ - **Christina Strelj**⁽²⁾ - **Martina Marchetti-Deschmann**⁽¹⁾

TU Wien, Institute of Chemical Technologies and Analytics, Vienna (1) - TU Wien, Atominstitut, Vienna (2)

Keywords: Multimodal Imaging, MALDI MS, μXRF , Sample Preparation, undecalcified tissue

Introduction

MALDI MSI visualizes distributions of a broad variety of biomolecules in tissue. μXRF is a non-destructive

elemental imaging method with high analyte specificity, and therefore a favourable choice for a multimodal imaging approach targeting elemental and molecular distributions. A methodology is presented allowing μ XRF and MALDI MSI analysis from the same undecalcified tissue section to visualize such distributions in tissue.

Methods

Chicken feet were purchased from a local market, digits were removed and embedded in a hydrogel containing NaCMC and gelatin. Cryosectioning of the phalanx was performed with a non-disposable tungsten carbide knife according to the Ullberg's method. An in-house built μ XRF spectrometer (Atominstitut, TU Wien) was used for XRF analysis and MALDI MSI experiments were carried out on a Waters Synapt G2 after applying matrix (dithranol) by sublimation.

Results

In this proof-of-principle study, we established and tested sample preparation strategies, allowing multimodal analysis of lipids (sphingomyelin and phosphatidylcholines) and elements relevant for bone structures as calcium, phosphorous and sulphur in the very same sample section of a chicken phalanx without tissue decalcification. The results of the investigation of such parameters as adhesive tapes supporting tissue sections, and sequence of the imaging experiments are presented. We show specific lipid distributions in skin, cartilage, muscle, nail, and intact morphology of bone by calcium and phosphorous images. Combination of molecular and elemental imaging was achieved, thus, providing now for the first time the possibility to gather MALDI MSI and μ XRF information from the very same sample without any washing steps omitting therefore analytical artifacts that inevitably occur in approaches using consecutive tissue sections.

Conclusions

Conversely to most previous multimodal (multisensor) investigations, our imaging experiments are carried out on the very same samples without any additional steps, thus, avoiding analytical artefacts. We showed that we are able to cut difficult samples properly and visualize lipid distributions together with biologically relevant elements like Ca, P or S. We performed imaging experiments with similar lateral resolution facilitating the correlation of the images and gaining morphological information that was not possible from either analytical method independently.

Novel Aspect

Combining MALDI MSI and μ XRF can benefit in researches where elemental and lipid interaction play an essential role.

WP-143 / COMPLEMENTARY BIOIMAGING TO DETERMINE NANOPARTICLE AND PHOSPHOLIPID DISTRIBUTION IN LUNG TISSUE

Ann-Christin Niehoff⁽¹⁾ - ***Dörthe Dietrich***⁽²⁾ - ***Michael Sperling***⁽²⁾ - ***Uwe Karst***⁽²⁾

Shimadzu Europa GmbH, Innovation Center, Duisburg (1) - University of Münster, Institute of Inorganic and Analytical Chemistry, Münster (2)

Keywords: Imaging, MALDI-MS, LA-ICP-MS, Multimodality, Nanoparticles

Introduction

The uptake of nanoparticles (NP) such as cerium oxide or zinc oxide NP via inhalation may pose a health hazard for humans. Therefore, knowledge about their distribution in the lung and their biological effects is needed. Here, we combined elemental and molecular imaging techniques to study the distribution of NP and changes of phospholipid distribution in the lung.

Methods

Different doses of nanoparticles were intratracheally instilled into rat lungs. The lungs were removed and cryo-sections were prepared. Laser ablation inductively coupled plasma mass spectrometry (LA-ICP-MS) experiments were performed to reveal the distribution of NP in the tissue. To analyze phospholipid distribution in positive and negative ion mode, matrix assisted laser desorption/ionization mass spectrometry experiments were carried out.

Results

To analyze the impact of nanoparticles on lung tissue, MALDI-MS experiments were performed. Different phosphatidylcholines were detected in the positive ion mode and verified via MS/MS experiments by characteristic fragment of m/z 184. In case of lung inflammation increasing signal intensities of PC(32:1) and PC(36:4) were detected. In the negative ion mode phosphatidylglycerols and phosphatidylinositols were detected and were confirmed via MS/MS experiments. Nanoparticle distribution patterns as well as the distribution of naturally occurring elements such as phosphorus were depicted by means of LA-ICP-MS. Accumulations of nanoparticles could be detected in alveolar areas of the lung.

Conclusions

The combination of elemental and molecular mass spectrometry by MALDI-MS and LA-ICP-MS shows the impact of nanoparticles on lung tissue. Depending on elemental composition, concentration, and period of incubation lung inflammation and fibrosis can occur. Alterations in phospholipid distribution were detected by MALDI-MS.

Novel Aspect

Results show that combining MALDI-MS and elemental imaging techniques can provide new insights into local changes caused by nanoparticles in tissues.

WP-144 / DISTINGUISHING THE CORE AND PENUMBRA BY LIPID PROFILES USING MASS SPECTROMETRY IMAGING IN MICE AFTER EXPERIMENTAL INFARCT INDUCTION

Nina Ogrinc Potocnik⁽¹⁾ - ***Inge Mulder***⁽²⁾ - ***Ludo Broos***⁽³⁾ - ***Arnoud Prop***⁽¹⁾ - ***Arn Maagdenberg***⁽²⁾ - ***Ron M. A. Heeren***⁽¹⁾

Maastricht MultiModal Molecular Imaging Institute (M4I), Maastricht University, Maastricht (1) - Leiden University Medical Center, Department of Neurology, Leiden (2) - Leiden University Medical Center, Department of Human Genetics, Leiden (3)

Keywords: Ischemic stroke, penumbra, MALDI-MSI, lipids, migraine

Introduction

Changes in lipidomic profiles have been linked to ischemic infarcts but little is known about these changes in the early phase of infarct evolution. Localization of species in the ischemic area in early stages is crucial for infarct development prognosis and possible clinical intervention. Mass Spectrometry Imaging was used to localize and identify specific lipid profiles at 3 different time points within the first 24 hours after experimental infarct (MCAO) in mice.

Methods

Four brains per time point at 4, 8 and 24 hours after MCAO were frozen using the in-situ funnel-freezing technique. DHB and Norharmane were sprayed at 15 mg/mL and 7 mg/mL in 2:1 (CHCl₃:MeOH), respectively with a TM-sprayer. In total over 120 sections were analysed with the Bruker Rapiflex MALDI TissueTyper™ at 50 μm spatial resolution. After MALDI-MSI sections were stained with a Nissl stain. High-resolution images were acquired at 24h post MCAO with nanoTOF II system.

Results

A linear transparency overlay, clustering pipeline and spatial segmentation revealed specific segments of the infarct area: the core, borderzone/penumbra (*i.e.* the area between the infarct core and healthy tissue) and healthy tissue. The borderzone became clearly visible at 8 hours post MCAO by $m/z = 965.5$ and 1045.5 that corresponded to phosphatidylinositol 4-phosphate (PIP1) and phosphatidylinositol 4, 5-bisphosphate (PIP2). Infarct evolution was clearly shown by 2D-renderings of multiple PC and Lyso-PC isoforms. Additionally, high-resolution secondary ion mass spectrometry (SIMS) was conducted to evaluate sodium/potassium ratios. These data revealed a significant increase in sodium species and a decrease in potassium distribution in the ischemic area compared to healthy tissue at 24 hours post MCAO. A more pronounced discrimination between the border zones and healthy regions was shown in a mouse model known for increased susceptibility to peri-infarct depolarizations and infarct vulnerability.

Conclusions

The observed changes deepen our understanding of disease mechanisms occurring in the first hours of stroke, and could lead to early detection of necrosis and could guide the development of new therapeutic targets to minimize infarct deficit.

Novel Aspect

Our findings may guide the development of early ischemic markers with the hope of a better prognosis concerning tissue fate in the penumbra.

WP-145 / THE MANAGEMENT OF HAEMOGLOBIN INTERFERENCE FOR THE MALDI-MSI ANALYSIS OF IN-VIVO THYROID BIOPSIES

Isabella Piga⁽¹⁾ - *Giulia Capitoli*⁽²⁾ - *Vanna Denti*⁽¹⁾ - *Silvia Tettamanti*⁽¹⁾ - *Stefania Galimberti*⁽²⁾ - *Fulvio Magni*⁽¹⁾ - *Fabio Pagni*⁽³⁾

University of Milano-Bicocca, Department of Medicine and Surgery - Clinical Proteomics and Metabolomics Unit, Veduggio al Lambro (1) - University of Milano-Bicocca, Department of Medicine and Surgery - Centre of Biostatistics for Clinical Epidemiology, Monza (2) - University of Milano-Bicocca, Department of Medicine and Surgery - Section of Pathology, Monza (3)

Keywords: Fine needle aspiration, thyroid lesion, haemoglobin interference, MALDI-MSI proteomics, thin prep.

Introduction

Fine Needle Aspiration biopsy (FNAB) is the gold standard exam to determine the malignant nature of thyroid nodules[1]. Contamination of FNAB samples with red blood cells is problematic for proteomics analysis, given that large amounts of haemoglobin (Hb) suppress other protein signals[2]. Hence, it is paramount to standardise the sample preparation of ex-vivo and in-vivo thyroid FNABs for proteomic MALDI-MSI analysis, in order to minimise Hb interference.

Methods

Human FNABs were collected and deposited onto conductive glass slides from both ex-vivo(n=9), surgically removed thyroid specimens, and in-vivo(n=12) thyroid specimens for intact proteins MALDI-MSI analysis. Three protocols were compared using ex-vivo biopsies collected from the same thyroid: a) conventional air-dried smear; b) cytological smear immediately fixed in ethanol; c) ThinPrep (TP) cytological preparation using a ThinPrep 2000 system.

Results

The spectral profiles of both ex-vivo and in-vivo conventional air-dried smears were characterized by high

inter-patient variability related to the abundance of the Hb peaks. In particular, the strong vascularization of some thyroid nodules is reflected in FNABs with a high content of Hb. The amount of Hb was markedly decreased in TP preparation with respect to both conventional air-dried and fixed smears. On the other hand, the absolute intensity of other protein signals, suppressed with the other two methods, were significantly increased in TP samples. Furthermore, the management of Hb interference of ex-vivo and in-vivo TP samples was comparable, indicating the opportunity to use in-vivo TP specimens for MALDI-MSI proteomic analysis and biomarker discovery. The MALDI-MSI approach combined with virtual microdissection permitted to extract specific protein signatures from different histotypes of both benign and malignant thyroid cell clusters.

Conclusions

The Thin Prep procedure for thyroid FNABs samples preparation combined with MALDI-MSI proteomic analysis allow us to obtain high-quality spectra, follicular cells specific protein profiles and to manage the haemoglobin interference. The application of this reproducible technique to in-vivo cytological samples can help cytopathologists in the diagnosis of thyroid nodules combining both morphological and proteomics information.

Novel Aspect

This study represents the first example of MALDI-MSI applied to ex-vivo and in-vivo thyroid FNABs, prepared using the ThinPrep preparation, for proteomic analysis.

References

1. Russ G, Bonnema SJ, Erdogan MF, Durante C, Ngu R, Leenhardt L. *European Thyroid Journal*, 6, 225-237 (2017).
2. Amann JM, Chaurand P, Gonzalez A, Mombley J, Massion PP, Carbone DP, Caprioli RM, *Clinical Cancer Research*, 12, 5142–5150 (2006).

Funding: This work was funded thanks to AIRC (Associazione Italiana per la Ricerca sul Cancro) MFAG GRANT 2016 - Id. 18445.

WP-146 / HIGH SPATIAL RESOLUTION IMAGING OF AGED HUMAN RETINA TISSUE

Kevin Schey⁽¹⁾ - ***David Anderson***⁽¹⁾ - ***Jeffrey Messinger***⁽²⁾ - ***N. Heath Patterson***⁽¹⁾ - ***Christine Curcio***⁽²⁾

Vanderbilt University, Biochemistry, Nashville (1) - University of Alabama-Birmingham, Ophthalmology, Birmingham (2)

Keywords: imaging mass spectrometry, retina, retinoids, lipids

Introduction

The human retina is made of morphologically distinct layers responsible for the absorption of photons and for sending signals necessary for vision to the brain. Imaging mass spectrometry (IMS) is a powerful tool for examining spatially disparate molecular signatures. We have applied IMS to human retina tissue to examine the spatial distribution of analytes in retina cell layers and to detect differences with age and disease.

Methods

Fixed and frozen human retina tissue from older donors (80-89 years) was sectioned to 12-14µm thickness and placed on ITO slides for IMS experiments. MALDI matrix was applied by sublimation and images were acquired using a modified Bruker Solarix 9.4T FTICR instrument. Data were analyzed with FlexImaging software. Registration between optical images and MALDI images was aided by autofluorescence microscopy and custom registration software.

Results

Fixed tissue sections provided superior morphology with reduced analyte signal intensity in contrast to frozen sections that had poorer morphology but provided higher signal intensity. Ion maps unique to many retina cell layers were obtained with 15 μm spatial resolution. Autofluorescence microscopy facilitated registration of IMS and microscopy via explicit mapping of IMS pixels to their MALDI laser ablation marks, allowing confident identification of retina cell layers in the IMS data. Molecular signatures were identified for the high density of Mueller cells within the fovea, for basal laminar deposits, and for the retinal pigment epithelium (RPE). The bis-retinoid metabolite, A2E, was detected in RPE cells of peripheral retina. A series of BMP lipids were detected throughout the RPE and photoreceptor outer nuclear layer. LC-MS/MS analysis is being used to identify novel compounds within the basal laminar deposits and within the fovea region. Unique signals were also observed within the inner and outer segments of photoreceptor cells.

Conclusions

IMS provides a powerful tool to examine molecular changes in distinct retinal cell layers with age and extracellular deposits in disease. Specific lipid and retinoid metabolites have been localized to the RPE and fovea regions of older human retinas.

Novel Aspect

Novel analytes identified in distinct retinal layers using high resolution IMS of human retina sections with characteristic aging features.

WP-147 / AMBIENT MASS SPECTROMETRY IMAGING ON WHOLE BODY SKIN FOR BIOMOLECULES

Jentaie Shiea

National Sun Yat-Sen University, Department of Chemistry, Kaohsiung

Keywords: skin, thermal desorption, ambient mass spectrometry, imaging, biomolecule

Introduction

Skin plays a vital role in protecting human being against the attack by pathogenic microorganisms. Biomolecules such neurotransmitters, amino acids and other metabolites are also released from skin through numerous glands underneath it. Although these biomolecules on skin are possibly potential disease biomarkers, current analytical techniques are incapable to efficiently identify them on skin, except blood drawing. In this study, thermal desorption-electrospray ionization/mass spectrometry (TD-ESI/MS), a novel ambient mass spectrometric technique[1], was used to characterize the skin compounds including neurotransmitters, amino acids, lipids and metabolites without doing any sample pretreatment. Since TD-ESI/MS owns the feature of rapidness, high-throughput sample analysis, molecular imaging of skin can be presented on specific chemical compounds. Herein, ambient mass spectrometry imaging was demonstrated for visualizing the distribution of biomolecules on whole body skin.

Methods

The TD-ESI/MS system consists of sampling probe, heated oven, electrospray device, and mass analyzer. A stainless steel probe was used to gently scrapes through the skin surface for ca 1 cm for sampling. It was then inserted into a heated oven to thermally desorb the analytes on it. The analytes were carried by a nitrogen stream down to a plume consisting of acidic methanol/water and ionized by reacting with the charged species (*e.g.* H^+ , H_3O^+ , etc.) in the ESI plume. Any residual sample on the probe was efficiently removed by burning it with a high-temperature gas torch for a few seconds. It took less than 30 seconds to complete an analysis.

Results

In typical case, approximate 2,000 samples were collected on skin within an hour and the analysis of these samples took less than 12 hours. The imagings were obtained using commercial available imaging software. Taking caffeine or medicine as an example, the analyte ion signal was mostly detected on the upper body especially forehead, fingertips and backside of the body. This may due to different sweating rates in different skin parts as well as the distribution of veins and capillaries over the body. Not only medicines and their metabolites but regular metabolites such as neurotransmitters, amino acids, fatty acids, cholesterol, diglycerides, triglycerides and squalene were also detected and their whole body skin imagings show that the distribution of metabolites on skin also depend on the structure and chemical property of the compounds.

Conclusions

Numerous biomolecules on the human skin was rapidly detected by in situ probe sampling and TD-ESI/MS analysis within a few seconds per sample. The ambient molecular imaging of squalene and cholesterol on whole body skin was also successfully established. In summary, we combined in situ probe sampling and TD-ESI/MS analysis, anon-invasive skin imaging technique for efficiently perform molecular imaging on the skin with high-throughput analysis was developed.

Novel Aspect

Thermal desorption-electrospray ionization/mass spectrometry was used to perform whole body skin ambient imaging on potential biomolecules.

References

1. Huang M. Z., Zhou C. C., Liu D. L., Jhang S. S., Cheng S. C., Shiea J., *Anal. Chem.*, 85, 8956-8963 (2013)

WP-148 / ALDOSTERONE SPECIFIC VISUALIZATION IN PRIMARY ALDOSTERONISM USING IMAGING MASS SPECTROMETRY

Shuichi Shimma⁽¹⁾ - *Emi Takeo*⁽¹⁾ - *Koshiro Nishimoto*⁽²⁾ - *Yuki Sugiura*⁽³⁾

Osaka University, Graduate School of Engineering, Osaka (1) - Saitama Medical University, Department of Uro-Oncology, Saitama (2) - Keio University School of Medicine, Department of Biochemistry, Tokyo (3)

Keywords: Imaging mass spectrometry, derivatization, structural isomer, steroid hormone

Introduction

Visualization of steroid hormones using MALDI-IMS is one of the active subjects for endocrine studies. Unfortunately, it is difficult to detect steroid hormones due to the low polarity. To overcome the problems, on-tissue chemical derivatization (OTCD) was developed [1], however the problem of structural isomers in steroid hormones is still remain. Herein, we introduce aldosterone specific visualization method using OTCD and tandem MS.

Methods

Frozen tissue was sectioned and thaw-mounted on conductive glass slides. Girard's reagent T (GirT) was used and the GirTsolution was sprayed on the tissue surface with an artistic air-brush. After one hour reaction at room temperature, α -CHCA was applied using two-step matrix application method. MALDI-IMS analysis was performed using LTQ XL linear ion trap mass spectrometer equipped (Thermo Scientific, Bremen, Germany).

Results

We first demonstrated discrimination of aldosterone and cortisone which is a structural isomer of aldosterone using authentic standard sample. Aldosterone and cortisone were successfully discriminated

by MS3 analysis. To confirm the utility of our method, next two experiments were performed in rat sodium diet model and human primary aldosteronism (PA), which is one of the causes of secondary hypertension. According to our experimental result, aldosterone and structural isomer of cortisone were clearly visualized separately. In PA patient tissues, we successfully visualized these hormones individually and compared with immunohistochemistry of aldosterone synthase. Our finding in this study was the expression level of aldosterone production enzyme and produced aldosterone were different in distribution.

Conclusions

Although some literatures described MALDI-IMS approach for steroid hormone imaging, structural isomers have not been fully considered. In this study, we have developed aldosterone specific visualization method. After validation of the method in rat adrenal gland, we successfully visualized aldosterone producing region in PA. We found that combination of different imaging results (MALDI-IMS and IHC) were valuable to elucidate diseases.

Novel Aspect

Combination of GirT and MS3 analysis allowed us aldosterone specific imaging in human adrenal gland, and it will provide more reliable information in diagnosis of PA.

References

1. D. F. Cobiceet al., *Anal Chem*, 85, 11576-11584 (2013).

WP-149 / MALDI-MS IMAGING IN THE SEARCH FOR PROTEOMIC INDICATORS OF RESPONSE TO THERAPY IN MEMBRANOUS NEPHROPATHY

Andrew James Smith⁽¹⁾ - ***Vincenzo L'Imperio***⁽²⁾ - ***Mariafrancesca Mazza***⁽¹⁾ - ***Elena Ajello***⁽²⁾ - ***Fabio Pagni***⁽²⁾ - ***Fulvio Magni***⁽¹⁾

Department of Medicine and Surgery, University of Milano-Bicocca, Clinical Proteomics and Metabolomics Unit, Veduggio al Lambro, Italy (1) - Department of Medicine and Surgery, University of Milano-Bicocca, Section of Pathology, Monza, Italy (2)

Keywords: clinical proteomics, high spatial resolution, membranous nephropathy, MS imaging, pathology

Introduction

Membranous Nephropathy (MN) is the most frequent cause of nephrotic syndrome in adults and the disease course is characterised by the “rule of third”, with one third of patients experiencing complete remission and the remaining experiencing relapses or progression [1]. Moreover, the therapeutic approach is not standardised, leading to further heterogeneity in terms of outcome [2], and a marker that can reliably predict the response of a patient a priori is lacking.

Methods

MALDI-MS imaging of renal biopsies taken from patients who differentially responded to immunosuppressive treatment (Ponticelli regimen) was performed. Initially, FFPE specimens (n=13) were analysed using an ultrafleXtreme™ MALDI-TOF-MS instrument, employing a pixel width and raster setting of 40 µm and 50 µm, respectively [3]. Analysis of additional specimens (n=7) was performed using a rapifleX MALDI TissueTyper™, employing a pixel width and raster of 10µm.

Results

Following Receiver Operative Characteristic (ROC) analysis, three signals were determined to have a statistical significance when comparing the two patient groups (*m/z* 1111, 1198 and 1303). The signal at *m/z* 1303 displayed the greatest discriminatory power and was observed to be of higher intensity in the glomeruli of the non-responding patients. The corresponding tryptic peptide was putatively identified as

macrophage migration inhibitory factor (MIF) and validated using immunohistochemistry, depositing on podocyte and parietal epithelial cells (PEC).

We also generated specific profiles of the mesangial cells and podocytes of the glomeruli whilst employing a 10µm pixel. This high spatial resolution imaging was then used to visualise the spatial localisation of those additional discriminatory signals detected (m/z 1111 and 1198) within the glomerular and tubular structures. These signals were putatively identified as tryptic peptides deriving from sonic hedgehog protein (SHH) and α -smooth muscle actin (SMA), respectively.

Conclusions

Despite much effort being made in recent years to understand the pathogenesis of MN, a biomarker able to predict the outcome of these patients following therapeutic treatment is still lacking. Here, we highlight three proteins that could differentiate between these MN patients (MIF, SHH and α -SMA) and thus represents a promising starting point in the search for protein markers to provide complimentary support in the routine prognostic assessment of MN patients.

Novel Aspect

This work represents the first study to employ high spatial resolution MALDI-MSI as a means of detecting tissue proteins with the potential to predict response to therapy in MN.

References

1. Ponticelli C. J. Nephrol., 20(3), 268–87 (2007)
2. Ponticelli C., Glasscock R. J., J. Am. Soc. Nephrol., 9, 609–16 (2014)
3. Smith A. et al. Biochim. Biophys. Acta, 1865, 865–874 (2017)

The research leading to these results has received funding from MIUR: FIRB 2007 (RBRN07BMCT_11), FAR 2014–2016 and in part by Fondazione Gigi & Pupa Ferrari Onlus.

WP-150 / PRECISION MEDICINE FOR DIABETIC NEPHROPATHY AND HYPERTENSIVE NEPHROSCLEROSIS: A MALDI-MSI PILOT STUDY

Andrew James Smith⁽¹⁾ - **Vadim Iablokov**⁽²⁾ - **Martina Stella**⁽¹⁾ - **Mariafrancesca Mazza**⁽¹⁾ - **Bram Heijs**⁽³⁾ - **Peter Van Veelen**⁽³⁾ - **Mariafrancesca Mazza**⁽¹⁾ - **Daniel Muruve**⁽²⁾ - **Fulvio Magni**⁽¹⁾

Department of Medicine and Surgery, University of Milano-Bicocca, Clinical Proteomics and Metabolomics Unit, Vedano al Lambro, Italy (1) - Department of Medicine, Cumming School of Medicine, University of Calgary, Calgary, Canada (2) - Center for Proteomics and Metabolomics, Leiden University Medical Center, Leiden, The Netherlands (3)

Keywords: chronic kidney disease, clinical proteomics, diabetic nephropathy, hypertension, MS imaging

Introduction

Diabetic nephropathy (DN) and hypertensive nephrosclerosis (HN) represent the most common causes of chronic kidney disease and many patients progress to end-stage renal disease [1]. Whilst both diseases are treated primarily through the management of hypertension, the latter group of patients have a more favourable outcome. Separating these two entities would allow for more appropriate resource allocation to those patients at highest risk of progression.

Methods

MALDI-MS imaging of renal biopsies taken from patients with DN and HTN was performed. These patients were also stratified into disease sub-groups on the basis of clinical parameters. FFPE specimens (n=18) were analysed using a rapiflex MALDI TissueTyper™, employing a pixel width and raster of 10µm and 20 µm,

respectively. Subsequent identification was performed using a Dionex UltiMate 3000 rapid separation (RS) LC nano system coupled with an Impact HDTM UHR-QqToF.

Results

Unsupervised Principal Component Analysis (PCA) highlighted that DN and HN could be clearly distinguished on the basis of their tryptic peptide profiles, with a cluster of signals having an altered intensity between the two disease groups. However, there was also some heterogeneity noted within each disease group. Upon further investigation, this heterogeneity was shown to be associated with the severity of each disease as defined by both clinical parameters and histological evaluation. Receiver Operative Characteristic (ROC) analysis highlighted a small panel of signals whose differential intensity was associated with the different stages of each disease. These signals were then correlated with a protein identification list obtained by complimentary nLC-ESI-MS/MS analysis, resulting in a putative cluster of proteins able to discriminate between the different stages of DN and HN.

Conclusions

A sensitive and specific protein marker that can reliably differentiate between diabetic and hypertensive causes of CKD would allow for targeted monitoring guidelines for clinicians. The cluster of putative proteins presented here represents a valuable starting point for a future study employing a larger cohort of patients and, if verified, would allow for more appropriate disease intervention and potentially improve the outcome of patients with DN and HN.

Novel Aspect

This complementary proteomics strategy may provide specific markers for diabetic and hypertensive causes of CKD and highlight those patients at highest risk of progression.

References

1. Jha V., Garcia-Garcia G., Iseki K., et al. *Lancet*, 382 (9888), 262-72 (2013)

WP-151 / HIGH-RESOLUTION AP-SMALDI MS IMAGING OF PARASITES

Bernhard Spengler⁽¹⁾ - **Patrik Kadesch**⁽¹⁾ - **Tobias Hollubarsch**⁽¹⁾ - **Anja Taubert**⁽²⁾ - **Carlos Hermosilla**⁽²⁾ - **Christoph Greveling**⁽²⁾ - **Thomas Quack**⁽²⁾ - **Stefanie Gerbig**⁽¹⁾ - **Kerstin Strupat**⁽³⁾ - **Mario Kompauer**⁽¹⁾ - **Max Müller**⁽¹⁾ - **Sven Heiles**⁽¹⁾

Justus Liebig University, Institute of Inorganic and Analytical Chemistry, Giessen (1) - Justus Liebig University, Institute of Parasitology, Giessen (2) - Thermo Fisher Scientific (Bremen) GmbH, Bremen (3)

Keywords: High-resolution MS imaging, neglected tropical diseases, parasites, 3D surface imaging, atmospheric pressure MALDI

Introduction

Biomolecular analysis by atmospheric-pressure (AP) MALDI MS imaging (MSI) at high lateral resolution (1) is a key issue for the investigation of small biological objects such as multicellular or protozoan parasites. Analysis of non-planar surfaces or tissue sections of parasite can be useful for, e.g., understanding invasion mechanisms. 3D-surface MSI2 allows describing the biochemistry, chemical communication and drug response of small objects.

Methods

An AP-SMALDI 5 AF ion source (TransMIT GmbH, Giessen, Germany) on a Q Exactive HF X orbital trapping mass spectrometer (Thermo Fisher Scientific, Bremen, Germany) was used to image two-dimensional planar sections and three-dimensional non-planar surfaces at a lateral resolution of 10 µm and below. High-speed acquisition up to 10 pixels per second was used at reduced mass resolution. Full pixel mode

enhanced signal intensities for low-abundant compounds.

Results

Schistosoma mansoni, a dioecious worm parasite causing schistosomiasis in humans and animals, was obtained from the Institute of Parasitology. 100 to 300 µm thick and approx. 1.5 mm long worms were investigated as couples or as individuals. Contact surfaces (the teguments) were differentially analyzed at the phospholipid and peptide levels, indicating molecular markers for males and females, and for their mating state.

Samples of apicomplexan protozoa in cell pellets and monolayers of bovine umbilical vein endothelial cells (BUVEC) were studied at the lipidome level to find markers of infection. More than 20 lipid biomarkers for infection were determined for each of the parasites *Besnoitia besnoiti*, *Eimeria bovis* and *Toxoplasma gondii*. Pixelwise autofocusing, provided by the 3D imaging ion source, was essential to reliably analyze the samples. Height adjustment in the micrometer range allowed obtaining topography-independent, concentration-related signal intensities from the surface of parasites and non-flat tissue sections.

Conclusions

High-resolution methods were developed for pinpointing biochemical features of parasites such as schistosomes or coccidia. Molecular communication between paired male and female schistosomes was characterized by 3D surface imaging of their contact zone, the tegument.

Semi-quantitative analysis was enabled by pixelwise autofocusing, resulting in topography-independent signal intensities.

Novel Aspect

Molecular 3D surface imaging of endogenous compounds of parasites at 5 to 10 µm lateral resolution

References

1. Kompauer M, Heiles S, Spengler B, Nature Methods 14, 90-96 (2017).
2. Kompauer M, Heiles S, Spengler B, Nature Methods 14, 1156-1158 (2017).

WP-152 / A MOLECULAR SNAPSHOT OF CLEAR CELL RENAL CELL CARCINOMA: MULTIMODAL MALDI-MSI AND NLC-MS/MS OF N-GLYCANS AND PROTEINS

Martina Stella⁽¹⁾ - ***Fanny Boyaval***⁽²⁾ - ***Andrew James Smith***⁽¹⁾ - ***Marco Grasso***⁽³⁾ - ***Giorgio Bovo***⁽⁴⁾ - ***Bram Heijs***⁽²⁾ - ***Fulvio Magni***⁽¹⁾

Univeristy of Milano Bicocca, Dept. of Medicine and Surgery, Proteomics and Metabolomics Unit, Monza (1) - Leiden University Medical Center, Center for Proteomics and Metabolomic, Leiden (2) - San Gerardo Hospital, Urology Unit, Monza (3) - Vimercate Hospital, Pathology Unit, Vimercate (4)

Keywords: Clear cell Renal Cell Carcinoma, MALDI-MSI, nLC-ESI-MS/MS, glyco-proteomics

Introduction

Renal cell carcinoma (RCC) is the most frequent form of kidney cancer and approximately 80% of cases are defined as clear cell RCC (ccRCC) [1]. Understanding molecular alterations associated with the onset and progression of the disease could provide a valuable starting point for the development of new target therapies. Therefore, combined high resolution mass spectrometric approaches were used to investigate the glyco-proteome of different ccRCC grades [2-3].

Methods

FFPE tissue samples from patients with a proven diagnosis of ccRCC were analysed by MALDI-MSI. A 9.4T MALDI-FTICR-MS platform was used to investigate the spatial distribution of proteolytic peptides, whereas

released *N*-glycans were analysed with a MALDI-TOF/TOF-MS platform from a consecutive tissue section. Finally, bottom-up proteomic using nLC-ESI-MS/MS, was performed on homogenized tissue to assign putative identifications to the molecules of interest.

Results

Considering the possibility of tumour grade-heterogeneity within a single lesion, traditional histology and molecular information obtained by MALDI-MSI were used to better define signals specific for each region of the tissue. MALDI-MSI analysis highlighted peptides signals able to discriminate between grade one and grade four lesions (e.g. m/z 944.71, m/z 1032.78, m/z 1325.99, m/z 1428.92). Nevertheless, preliminary data on the characterization of glycosylation offered an additional layer of information regarding the glyco-phenotype of different areas of the tissue. Finally, nLC-ESI-MS/MS analyses provided a list of putative identifications (protein and glycoproteins) for those signals that differed significantly among the four grades. In particular, a putative identification was assigned to the four aforementioned signals and were tryptic peptides of H2, H3, H4 and Vimentin, respectively. The identification of other signals of interest and their released glycans is currently under investigation.

Conclusions

Proteomics represents a highly specific approach in the investigation of molecular changes associated with tumour lesions. Furthermore, the study of aberrant glycosylation is bringing consistent advantages in the understanding of ccRCC progression. All those findings represent a first attempt to highlight molecular differences among ccRCC grades and could provide further insights into the molecular mechanisms implicated in ccRCC progression.

Novel Aspect

The combination of high resolution MALDI-MSI and nLC-ESI-MS/MS is a highly effective approach for studying the proteome and the glycome of different tumour grades.

References

1. B. Ljungberg, K. Bensalah, S. Canfield, S. Dabestani, F. Hofmann, M. Hora, M. A. Kuczyk, T. Lam, L. Marconi, A. S. Merseburger, P. Mulders, T. Powles, M. Staehler, A. Volpe, A. Bex, Eur. Urol. 67, 913–924(2015)
2. Pinho SS, Reis CA, Nat Rev Cancer. 15, 540-55 (2015)
3. J. Dagher, B. Delahunt, N. Rioux-Leclercq, L. Egevad, J. R. Srigley, G. Coughlin, N. Duglison, T. Gianduzzo, B. Kua, G. Malone, B. Martin, J. Preston, M. Pokorny, S. Wood, J. Yaxley, H. Samaratinga, Histopathology. 71 (2017)

Acknowledgments

The research leading to these results has received funding from the MIUR: FIRB 2007 (RBRN07BMCT_11), FAR 2014–2016; and in part by Fondazione Gigi & Pupa Ferrari Onlus.

WP-153 / EXPLORING METABOLIC DEREGULATION IN MOUSE MODELS OF BREAST CANCER WITH MULTIMODAL MASS SPECTROMETRY IMAGING

Adam Taylor⁽¹⁾ - ***Tingting Fu***⁽¹⁾ - ***Bin Yan***⁽¹⁾ - ***Alex Dexter***⁽¹⁾ - ***Teresa Murta***⁽¹⁾ - ***Spencer A. Thomas***⁽¹⁾ - ***Yulia Panina***⁽²⁾ - ***Ala Al-Afeef***⁽¹⁾ - ***Efstathios Elia***⁽¹⁾ - ***Chelsea J. Nikula***⁽¹⁾ - ***Robin Philip***⁽¹⁾ - ***Rory T. Steven***⁽¹⁾ - ***Simon Barry***⁽³⁾ - ***Kevin Brindle***⁽⁴⁾ - ***Ian S. Gilmore***⁽¹⁾ - ***Richard J. A. Goodwin***⁽⁵⁾ - ***John Marshall***⁽⁶⁾ - ***George Pouligiannis***⁽⁷⁾ - ***Owen J. Sansom***⁽⁸⁾ - ***Zoltán Takáts***⁽⁹⁾ - ***Mariia O. Yuneva***⁽²⁾ - ***Josephine Bunch***⁽¹⁰⁾

National Physical Laboratory, National Centre of Excellence in Mass Spectrometry Imaging, Teddington (1) - The Francis Crick Institute, Oncogenes and Tumour Metabolism Laboratory, London (2) - AstraZeneca, Oncology iMED, Cambridge (3) - University of Cambridge, CRUK Cambridge Institute, Cambridge (4) -

AstraZeneca, Drug Safety and Metabolism, Cambridge (5) - Queen Mary University of London, Barts Cancer Institute, London (6) - Institute of Cancer Research, Division of Cancer Biology, London (7) - CRUK Beatson Institute, Colorectal Cancer and Wnt Signalling, Glasgow (8) - Imperial College London, Department of Surgery and Cancer, London (9) - National Physical Laboratory, National Centre of Excellence in Mass Spectrometry Imaging, London (10)

Keywords: MALDI, DESI, imaging, cancer, metabolism

Introduction

Overexpression of Myc and ErbB2 are frequent inducers of human breast cancers with distinct metabolic profiles [1]. Such metabolic deregulation is one of the hallmarks of cancer [2]. Mass spectrometry imaging (MSI) provides the potential to understand metabolic changes between genetic phenotypes, and within tumor regions; to expose therapeutic vulnerabilities. We performed a multimodal MSI study of genetically engineered mouse models of breast cancer.

Methods

Mouse normal mammary glands, Myc- and ErbB2-induced mammary gland tumours were snap-frozen and cryosectioned at 10 µm. MALDI MSI was performed on a Waters Synapt G2-Si at 45 µm pixel size using 9-AA or DHB matrix. DESI MSI was performed on a Waters Xevo G2-XS at 75 µm pixel size with 95:5 MeOH:water as the spray solvent. Data were analyzed using PCA, NMF, k-means clustering and deep learning t-SNE. Database matching was performed to tentatively assign ions of interest.

Results

Commonly applied analysis techniques PCA, NMF and k-means clustering proved complementary in revealing heterogeneity between subtypes and within individual tumours in an unsupervised manner. A novel deep learning t-SNE method enables visualisation of diverse chemical environments in a dataset that would be prohibitively large for traditional t-SNE methods. Consideration of spectral factors and loadings from multivariate analysis led to the identification of peaks that delineate chemical heterogeneity. Database matching was used for tentative identifications of molecules including lipids, fatty acids and small metabolites. Information provided by multivariate segmentation was broadly consistent across techniques. MALDI and DESI were found to be highly complementary, with a number of molecules uniquely identified by each technique alone. Future work will integrate MALDI and DESI MSI into a data guided approach for high spatial and mass resolution SIMS imaging using 3D OrbiSIMS.

Conclusions

We show the power of a multimodal MSI pipeline to delineate and define molecular heterogeneity within and between Myc- and ErbB2-induced models of breast cancer. MALDI and DESI MSI are shown to be complementary in their coverage, demonstrating the importance of multimodal MSI. This work provides a blueprint for the design and application of further multimodal and multiscale MSI studies.

Novel Aspect

A multimodal, multiscale MSI pipeline for the study of metabolic deregulation in genetically engineered mouse models of breast cancer is described.

References

1. Nair R. et al., *Oncogene*, 33, 3992-4002 (2014)
2. Hanahan D., Weinberg R. A., *Cell*, 144, 646-674 (2011)

WP-154 / NEW FRONTIERS OF MALDI-IMS AND LIPIDOMICS: AN UNTARGETED APPROACH TO SKIN LIPIDS MAPPING

Sara Tortorella⁽¹⁾ - ***Fang Xie***⁽²⁾ - ***Paolo Tiberi***⁽³⁾ - ***Gabriele Cruciani***⁽⁴⁾ - ***Stephen Castellino***⁽²⁾

Molecular Horizon, R&D, Bettona, PG (1) - GlaxoSmithKline, Bioimaging, In Vitro In Vivo Translation UM, King of Prussia, PA (2) - Molecular Discovery, R&D, Borehamwood (3) - University of Perugia, Department of Chemistry, Biology and Biotechnology, Perugia (4)

Keywords: MSI, lipidomics, chemometrics, skin, bioimaging

Introduction

Skin is large and integrated metabolic organ highly specialized in the synthesis of lipids, and skin untargeted lipidomics aim to facilitate the understanding of possible mechanisms of diseases with cutaneous pathological phenotype or pharmacokinetics of topically administered drugs [1,2].

Skin lipids present a specific distribution that might be exploited for the structural characterization of skin, thus obtaining a mapping that would guide any future detection of anomalies in the lipid profile. To this end, development of innovative high resolution and high-sensitivity analytical approaches characterized by speed, robustness, and high metabolite coverage is essential to empower the information associated with the skin lipidome. In this context, in addition to the powerfulness demonstrated by mass spectrometry in “omics” studies, mass spectrometry imaging (MSI), a molecular analytical technology suited for simultaneously measuring multiple analytes (including lipids, drugs and their metabolites) directly from intact tissue sections, adds a new dimension and holds potential for bringing lipidomics to the spatial molecular characterization level [3].

However, such an aim represents a remarkable computational challenge due to the size and analytical complexity of a typical MALDI MSI lipidomic analysis that, especially when high spatial resolution is used to characterize the tissues, increases the problem to the scale of Big Data. Therefore, a smart combination of bioinformatics and chemometric tools for image analysis and data mining is needed to effectively support MALDI MSI skin lipidomics, but adequate and user-friendly solutions comprising all of the steps required for this aim (from peak picking, to lipid identification via multivariate statistical analysis) are still lacking.

Here we propose a novel and semi-automatic workflow to map skin lipids that, starting from multivariate statistics, leads to skin lipids mapping and biomarkers discovery in an untargeted manner under orientation of histological structures by MALDI MSI.

To prove the efficacy of the proposed workflow, both minipig, a representative preclinical species for skin studies, and human skin samples were analysed with unprecedented spatial resolution (10 μm and 25 μm) under positive as well as negative ionization mode. Data were then processed by means of multivariate statistical analysis with the aim of identifying differences among skin features.

As a result, a relatively comprehensive set of specific lipids markers for epidermis, hair follicles, and for the first time, for sebaceous glands is identified. Structure identification of identified markers is confirmed by post-MALDI MSI H&E stained optical images and validated by additional samples.

Methods

Both minipig and human skin samples were collected and stored in -80°C before sectioning. The samples were cryosectioned at a thickness of 6 μm (minipig) or 10 μm (human) by a Leica Cryostat (Leica Biosystems Inc., Buffalo Grove, IL, USA) at -20 to -25°C . Sections for MALDI MSI were applied of matrix (DHB for positive and 9-AA for negative ionization mode) by a TM sprayer (HTX Technologies, Carrboro, NC, USA) upon optimization of spray conditions. MALDI MSI was performed using a Solarix 7T Fourier transform ion cyclotron resonance mass spectrometer (BrukerDaltonics, Billerica, MA).

Images were acquired at spatial resolutions of 10 μm and 25 μm in positive and negative ionization mode. Mass spectra were acquired in full scan mode (m/z 200-1500). All ion images were generated using FlexImaging v 4.0 software (BrukerDaltonics, Billerica, MA) from the raw data. Regions of interest (ROIs) were defined under the guidance of hematoxylin and eosin (H&E) staining as: epidermis, EP, papillary dermis, PD, reticular dermis, RD, hair follicle, HF, sebaceous glands, SG, sweat glands, SWG. At least three different skin portions were chosen for each ROIs. Average spectra from the different ROIs were exported from FlexImaging v 4.0 software. LipoStar software [4] (Molecular Discovery, London, UK), adapted for handling MALDI MSI datasets, was used for peak detection, samples alignment as data matrix (rows: samples or ROIs, columns: m/z values), supervised (partial least squares discriminant analysis, PLS-DA) and

unsupervised (principal component analysis, PCA) multivariate statistical analysis and lipid identification.

Results

The skin has three main layers: the epidermis, dermis, and hypodermis. Epidermis (EP) is the surface layer in contact with the exterior, while the dermis (D) is divided into two further substructures: papillary dermis (PD), immediately adjacent to the EP, and the reticular dermis (RD). Additional skin appendages are present, namely: Hair follicles (HF), sebaceous glands (SG), and sweat glands (SWG). In order to investigate differences among substructures, tissues were analysed at high spatial resolutions (10 μm and 25 μm) under positive as well as negative ionization mode.

Data were pre-processed as described in methods, and the obtained aligned data matrix underwent both supervised and unsupervised multivariate statistical analysis.

PCA analysis shows that groups of lipid scan be identified to distinguish the different layers and structures at a spatial resolution of 25 μm for both minipig and human skin samples. For instance, a lipid was identified to be specifically localized in the epidermis and the external root sheath of hair follicles. Considering that the external root sheath consists of epithelial cells and finally merges with epidermis, this identification is consistent with histological interpretation. Another lipid was found to be specifically expressed in hair follicles that contain a shaft.

Moreover, in human skin samples, SWG and SG substructures were characterized, and by means of PLS-DA it was possible to identify lipids specifically located on those substructures. Noteworthy, a list of lipids specifically located in SG was identified for the first time.

At the higher resolution of 10 μm , it was possible to identify lipids discriminating EP from PD and substructures of SG in minipig and human skin.

Lipids specific for each substructure were tentatively annotated with LipoStar software and post-MALDI MSI H&E stained optical images were collected to confirm the structure identification. When isobars were identified, an additional inquiry to a database of lipids previously identified in skin samples allowed for a priority match assignment.

The obtained list of skin substructure-specific lipids provided a lipid mapping that will be a valuable help in future skin lipidomic investigations.

Conclusions

In this work we presented a novel high-throughput cheminformatic workflow to efficiently map skin lipids and ultimately achieve untargeted lipidomics with MALDI MSI.

With the aim of mapping skin lipids, both minipig and human skin samples were analyzed at 10 μm and 25 μm spatial resolution in positive as well as negative ionization mode. These two spatial resolutions allowed distinguishing among substructures of different size. Starting from the selection of different ROIs according to H&E staining, multivariate statistical algorithms (PCA, PLS-DA) were applied to enlightening analogies/differences among the different skin substructures (stratum corneum, epidermis, dermis, hair follicles, and for the first time, sebaceous glands). In the present work we focused our attention on finding substructure-specific lipids, that is, lipids that exclusively belong to a specific substructure. A first annotation of those lipids was also provided, resolving isobars conflicts by means of a database of lipids previously identified in skin.

Big challenges still remain, and these should guide the design of more sophisticated and fully integrated (from peak detection to pathway analysis, via lipid identification) workflow.

However, by designing a workflow to successfully bringing lipidomic studies to spatial molecular characterization, we believe that these results will pave the way to establish a data-driven histology interpretation of MALDI MSI analyses and ultimately achieve functional elucidation to support drug and biomarkers discovery.

Novel Aspect

With the aim of performing for the first time an untargeted lipidomic analysis of skin by means of MALDI MSI, a novel cheminformatic workflow was designed to fish for lipids that are specifically localized in given skin substructures. Both minipig and human skin samples were analysed by MALDI MSI at unprecedented

spatial resolution (10 μm and 25 μm) in positive and negative ionization mode. As a result, a relatively comprehensive set of specific lipids markers for epidermis, papillary dermis, hair follicles, and for the first time, for sebaceous glands is identified. To the best of our knowledge, this is the first time that histology interpretation of MALDI MSI analyses is achieved in a completely data-driven fashion. These results will guide further pathway analyses to achieve unprecedented functional interpretation, to comprehensively map skin lipids and ultimately support drug and biomarker discovery.

References

1. Lee D. D., Zavadil J., Tomic-Canic M., Blumenberg M., Epidermal Cells: methods and Protocols (ed. Turksen, K.) 193–223 (Humana Press) (2010).
2. Camera E., Ludovici M., Tortorella S., Sinaglia J. L., Capitanio B., Goracci L., Picardo M., Journal of Lipid Research 57, 1051-1058 (2016).
3. Schwamborn K., Caprioli R. M., Molecular Oncology 4, 529-538 (2010).
4. Goracci L., Tortorella S., Tiberi P., Pellegrino R. M., Di Veroli A., Valeri A., Cruciani G., Analytical Chemistry 89, 6257-6264 (2017).

WP-155 / MALDI-MASS SPECTROMETRIC IMAGING FOR THE INVESTIGATION OF METABOLITES IN LINUM USITATISSIMUM L. ROOTS INOCULATED WITH FUSARIUM OXYSPORUM

Gleb Vladimirov⁽¹⁾ - *Nataliya Melnikova*⁽²⁾ - *Alexey Dmitriev*⁽²⁾ - *Nadezhda Bolsheva*⁽²⁾ - *Roman Novakovskiy*⁽²⁾ - *Alexander Zhrebker*⁽¹⁾ - *Artur Yablokov*⁽³⁾ - *Evgeny Nikolaev*⁽¹⁾

Skolkovo Institute of Science and Technology, Skoltech Center for Translational Biomedicine, Moscow (1) - Engelhardt Institute of Molecular Biology RAS, Engelhardt Institute of Molecular Biology RAS, Moscow (2) - Institute of energy problems of chemical physics RAS, Moscow (3)

Keywords: Linum usitatissimum, Metabolites, MALDI, Orbitrap, Imaging

Introduction

Flax (*Linum usitatissimum* L.) is used for fiber, oil, and nutraceutical production. *Fusarium oxysporum* f. sp. lini is the pathogen that causes fusarium wilt in flax. Epidemics of the disease can result up to 100% loss in yield. Some of flax genotypes are resistant to *F. oxysporum*, but the mechanisms of resistance are still unclear. Comparative analysis of metabolic changes in flax plants of resistant and susceptible cultivars in response to the pathogen infection is essential for understanding the role of metabolites in resistance and can be useful for development of antifungals.

Methods

Seven-day-old seedlings of one susceptible (AP5) and two resistant (Dakota and #3896) flax cultivars were inoculated with spore suspension of *F. oxysporum* pathogenic isolate #39 from the phytopathogen collection of the All-Russian Research Institute for Flax. For control plants, sterile water was used instead of *F. oxysporum* spore suspension. After 48 hours from inoculation, root tips were collected. After that, roots were embedded in gelatin and gently frozen on dry ice. The frozen tissue was then sectioned into 25- μm slices using a cryostat MK-25M (adapted for work with blades Leica DB80LS) at -20°C . The sections were thaw-mounted onto ITO slides. Matrix (40 mg/mL DHB in 50:50 water:acetonitrile) was applied using an airbrush (iwata micron cm-b2). A MALDI-Orbitrap mass spectrometer (Thermo Scientific Q-Exactive orbitrap with MALDI/ESI Injector from Spectrograph, LLC) that was equipped with a 355 nm laser (spot diameter of 20 μm) was used in positive ion mode for imaging and MS/MS [1]. Spectra were obtained in mass range of m/z 100–1000, a mass resolution of 140,000. The tissue region to be imaged and the raster step size were controlled using the Spectrograph MALDI Injector Software. To generate images, the spectra were collected at 35- μm intervals in both the x and y dimensions across the surface of the sample. Ion images were generated from raw files (obtained from orbitrap tune software) and coordinate files (obtained from MALDI

Injector Software) by image insight software from Spectrograph LLC. MS/MS collision-induced-dissociation (CID) fragments were collected by isolating each m/z of interest and manually adjusting the collision energy for each compound. The “metabolites of interest” were identified by searching the accurate mass obtained with the Orbitrap instruments and the MS/MS data using MetFrag.

Results

Metabolites in flax roots of susceptible and resistant cultivars inoculated with *F. oxysporum* or under control conditions were identified. For evaluation of metabolic responses to the fungal infection, comparative analysis was performed, and metabolites that were induced in resistant and susceptible flax genotypes after inoculation with *F. oxysporum* were identified. Resistance-related metabolites were determined. In the present study, we used susceptible and resistant flax cultivars, for which transcriptome analysis under the same conditions (inoculation with *F. oxysporum* isolate #39 and control) was performed by us previously, and diverse patterns of differentially expressed genes were revealed.

Conclusions

Complex analysis of transcriptome and metabolome data in flax plants with distinct resistance to fusarium wilt can be a promising tool for determination the mechanisms of resistance to *F. oxysporum*.

Novel Aspect

Metabolite profiling was performed for inoculated with *Fusarium oxysporum* flax roots of susceptible and resistant to the pathogen cultivars.

References

1. Belov M.E., Mikhail E., et al. "Design and Performance of a Novel Interface for Combined Matrix-Assisted Laser Desorption Ionization at Elevated Pressure and Electrospray Ionization with Orbitrap Mass Spectrometry" *Analytical Chemistry* 89, 7493-7501 (2017).

WP-156 / IDENTIFICATION OF ASBESTOS FIBRES IN A MALIGNANT MESOTHELIOMA MODEL

Oana Voloaca⁽¹⁾ - ***Laura Cole***⁽¹⁾ - ***Malcolm. R. Clench***⁽¹⁾ - ***Amy Managh***⁽²⁾ - ***C. Greenhalgh***⁽²⁾ - ***Sarah Haywood-Small***⁽¹⁾

Sheffield Hallam University, Biomolecular Science Research Centre, Sheffield (1) - Loughborough University, Department of Chemistry, Loughborough (2)

Keywords: Mesothelioma, LA-ICP-MS, Mineral Fibres, Imaging

Introduction

Malignant Mesothelioma (MM) is an aggressive cancer of the mesothelium associated with exposure to specific Mineral Fibres (MF) [1]. There is an urgent need to develop methods to clearly identify and quantify the MF within biological samples [2].

The current project aim was to develop a biologically relevant MM model and use LA-ICP-MSI to determine the element distribution of MF within samples. Development of this technique would untimely lead to early MM diagnosis.

Methods

MM models were developed using immortalised cells (NCIH28, MSTO211H, NCIH1975, Met5A) cultured in both 2D monolayers and 3D spheroids. Cells were exposed to various MF (concentration 0-25 $\mu\text{g/mL}$). Viability of the cultures was measured using Alamar Blue[®] and Hoescht/Propidium Iodide staining. MALDI and LA-ICP-MS images were obtained at post-exposure following cell cytospin and matrix coating. The spheroids were cryosectioned and fixed onto slides prior to the imaging.

Results

The MM models were successfully developed in both 2D monolayers spiked with MF and 3D spheroids with MF embedded within the cells. The Hoescht/Propidium Iodide staining showed an apoptotic core to the spheroids which translates to the in vivo tumour environments.

The cell viability results after MF treatment declined with exposure period, MF dose and most importantly MF type, which is in line with published literature [3]. Cell viability data of 2D models were also supported by Hoescht/Propidium Iodide images.

Some preliminary data has been acquired using LA-ICP-MS imaging, proving that this analytical method has high potential in identifying the MF within the cell cytopins and can become a useful tool for clinical diagnosis of MM. Initial data has also shown that sections of the 3D MM model can be analysed by LA-ICP-MS and MALDI MSI to determine spatial distribution of lipids and proteins surrounding MF.

Conclusions

LA-ICP-MS has recently emerged as one of the most versatile imaging tools in biomedical research and the current study has showed great potential for its future integration in clinical settings. Method optimisation is needed before moving on to more complex 3D structures and histopathological lung tissue biopsies. MALDI MSI and flow cytometry will also be employed in future analysis of protein expression within MM samples and to support current data.

Novel Aspect

This project is the first attempt to develop a diagnosis method that uses LA-ICP-MSI to identify the full arrangement of MF at the cellular level within MM samples.

References

1. Carbone C., Yang H., Clinical cancer research: an official journal of the American Association for Cancer Research 18, 598-604, 2012.
2. Chernova et al., Cell Death and Differentiation, 23, 1152–1164, 2016.
3. Hodgson JT., Darnton A., The Annals of occupational hygiene, 44, 565-601, 2000.

WP-157 / TWO DIMENSIONAL METABOLIC PROFILE OF BILE ACIDS AND LIPIDS IN FIBROSIS OF NASH-MODEL LIVER USING MALDI IMAGING AT 20 MM SPATIAL RESOLUTION

Yuzo Yamazaki⁽¹⁾ - ***Shinichi Yamaguchi***⁽²⁾ - ***Tomonori Oshikawa***⁽³⁾ - ***Daisuke Miura***⁽⁴⁾ – ***Makoto Yamazaki***⁽⁵⁾

Shimadzu Corp., Global application development center, Kyoto (1) - Shimadzu Corp., MS business unit, Kyoto (2) - Shimadzu Corp., Research & Development Department, Kyoto (3) - Kyusyu University, Innovation Center for Medical Redox Navigation, Fukuoka (4) - Mitsubishi Tanabe Pharma Corporation, DMPK Research Laboratories, Toda (5)

Keywords: MALDI imaging, NASH, lipid, bile acid

Introduction

Non-alcoholic steatohepatitis (NASH) is an advanced form of Non-alcoholic fatty liver disease (NAFLD), and is recognized to often progress to cirrhosis and hepatocellular carcinoma[1]. A liver excised from an animal model raised to forcibly express NASH is quite useful to study the disease, because pathological findings of the organ are similar to the one from human patients. Whereas some metabolomics methods using LC- or GC-MS have been applied to the analysis of the organ, visualization of metabolites on pathological section is necessary to discuss disease states with relevance to various staining. We will present a metabolic profile of biochemical compounds in some pathological states, which is possibly related to a dynamics of the fibrosis in liver, by using MALDI imaging.

Methods

Liver tissues were removed from 4w-, 8w-, and 16w-old male of NASH-model rat induced by high fat cholesterol diet. Dried liver sections were subject to MS imaging and some pathological staining (hematoxylin-Eosin, α SMA, sirius red, azan). Three sections with three corresponding controls were mounted onto an ITO-coated slide glass. 9-aminoacridine was sublimated onto the tissues at 0.5 μ m thickness with a commercial instrument (iMLayer, Shimadzu Corp., Japan). After sublimation, a recrystallization of the coated matrix was conducted under vapor of a solvent consisted with methanol and water. MS imaging at 20 μ m spatial resolution were performed with a MALDI-TOFMS (MALDI-7090, Shimadzu Corp., Japan) in negative mode at 2 kHz laser repetition, after a diameter of ablation by laser was validated. Newly developed in-house softwares were applied to visualization and analysis of all data.

Results

Three liver tissues, 4w, 8w, 16w-treated, and each corresponding control tissues were subject to the high through-put MS imaging simultaneously, the advantage of which is that distributions of biochemical compounds can be recognized undoubtedly under the same matrix coating and laser irradiation. Targeted m/z of biochemical compounds obtained from LC- or GC-MS was visualized with in-house software. As typical case, since an increase of PS(18:0,18:1) and PI(16:0,18:2) in 8w- and 16w-treated liver was indicated with LC-/GCMS, m/z of both lipids were visualized to be compared with stained sections. Whereas the PS was extremely localized at fibrosis area, which was stained with azan and sirius red deeply, the PI was found at nodules in the liver tissues. Bile acids and their derivatives were also chosen as target molecules and visualized. Interestingly, whereas taurocholate and taurochenodeoxycholate distributed at nodule areas widely in 8w- and 16w-treated liver, there were significant narrow areas in the nodule, where only the taurocholate or the taurochenodeoxycholate were accumulated specifically. On the other hand, glycocholate were localized in fibrosis areas that were increasing along with the treated week. Two-dimensional operation included in the newly developed software was applied to find more precise relevance between a progress of fibrosis and a distribution of bile acids. Calculated images of NASH liver tissues obtained by dividing m/z 405 (glycocholate) by summation of m/z 498 and 514 (taurochenodeoxycholate and taurocholate respectively) indicated that a distribution of glycocholate was correlated with fibrosis areas stained with azan in detail. In general, fibrosis in liver is attributable to excessive production of collagen and extracellular matrix in hepatic stellate cell after liver injury by stress. Specific 2D profiles of bile acids and lipids shown here could be caused with a change of global metabolism by the stress. A production of actual NASH liver in lab animals is still matter due to long necessary months. Raising the animals with high fat cholesterol diet is useful countermeasure to obtain almost the same pathology as the one in actual disease, however, pathological findings have been only the way to examine both tissues[2]. The 2D profiles obtained here is applicable to validate NASH-model at molecular level.

Conclusions

Differently treated NASH-model livers were studied using MALDI imaging to discuss a relevance with a pathological progress. Not only a phospholipid but also both primary and secondary bile acids were mapped in fibrosis area that increases with the HFC-treated week.

2D profiles of those molecules targeted with metabolomics analysis indicate a validation of NASH-model tissue could be possible with MALDI imaging, which might contribute to investigation and development of novel medicine for NASH.

Novel Aspect

MALDI imaging was applied to obtain metabolic profiles of bile acids and lipids in fibrosis of NASH-model liver.

References

1. Cohen JC., Science, 332, 1519-23(2011).
2. Horai, Y., et al., Int J Exp Pathol., 97, 75-85(2016).

WP-158 / MASS SPECTROMETRY IMAGING FOR DIRECT VISUALIZATION OF ENVIRONMENTAL STRESS INDUCED PHYTOHORMONES *IN-SITU*

Chao Zhang - Karel Doležal - Ondřej Novák

Laboratory of Growth Regulators, Palacký University Olomouc, Olomouc

Keywords: Plant wound stress, Jasmonic acid, abscisic acid, MALDI-MS Imaging

Introduction

Phytohormones are low concentration small molecules, but necessary for biological activities and development of plants, especially in metabolite regulations responding to the environmental stress. Previously, abscisic acid (ABA), Jasmonic acid (JAs) and other stress-induced phytohormones have been identified and quantified from treated leaves, but their spatial distributions on leaves are still unknown (Verma et al., 2016). Mass spectrometry imaging (MSI) is an advanced technique for small molecule studies which has been successfully applied in plant research. This research aims to highlight unique phytohormones patterns in-situ induced by different environmental stress and their relationships with plant developments. The visualization of phytohormones via DESI/MALDI MSI will not only complement to target profiling results achieved in UHPLC–MS/MS analysis but also deliver a better understanding on the stress-induced signaling of phytohormones

Materials and Methods

Arabidopsis thaliana plants harvested at the 24 day with half of the fully developed rosette leaves and treated on the one side of the central vein to induce environmental stress. Local stressed leaves were collected after 24h and removed cuticle with chloroform in order to visualize targeted phytohormones in-situ, distal leaves were used as the untreated control. DESI-IMS was performed using a SYNAPT G2-Si Q-TOF instrument equipped with a DESI ion source (Waters). The spray solvent (90% methanol, v/v) was introduced with a flow rate of 3 μ L/min using a Harvard syringe pump. A voltage of 3kV and 0.5MPa nitrogen was performed on the DESI sprayer. Leaf samples were scanned with 300 μ m resolution at 300 μ m/s, and collecting the mass spectrum every 0.5s. Signals at m/z 50–1200 were detected in the negative ion sensitive mode and ion images were constructed using the HDImaging software (Waters). Internal standards were performed using 1pmol/ μ L ABA and JA standards in this measurement. The ABA or OPDA standards were diluted to 1pmol/ μ L in 100% methanol and spot 1 μ L onto a glass slide. The MS/MS spectrum of each compound was measured using a SYNAPT G2-Si instrument equipped with the ion source (Waters). The spray solvent (90% methanol, v/v) introduced with a flow rate of 3 μ L/min using a Harvard syringe pump. A voltage of 3kV and 0.5MPa nitrogen was performed on the DESI sprayer. The optimized CE voltages for the fragmentations were 12eV for ABA and 18eV for JA. The chemical structures of the major fragment ions were compared with records in the online database.

Results

Intact masses of JA, ABA and other phytohormones were detected in the DESI MS spectra. 2. Spatial distributions of phytohormones on leaf samples were highlighted in their ion intensity maps. 3. DESI-MS imaging results showed unique hormonomics patterns in different wound healing stages. 4. PCA analysis and hypothesis tests suggested levels of JA and ABA are correlated with other phytohormones.

Conclusions

* DESI mass spectrometry imaging could identify and matched phytohormones profiles acquired in UHPLC–MS/MS from their intact positions.* The distribution differences of phytohormones near treated regions were highlighted in DESI-MSI.* Mass spectrometry imaging can be used as a discovery tool which could reveal new stress-related hormones in plant.* ABA and JA were found strongly correlated with plant metabolism alternations under different environmental stress.

WP-159 / A TOF-SIMS STUDY OF LIPID CHANGES IN E. COLI MUTANTS WITH IMPAIRED PLASMID TRANSFER CAPABILITY

Kelly Dimovska Nilsson - Martin Palm - Anne Farewell - John S. Fletcher

University of Gothenburg, Department of Chemistry and Molecular Biology, Gothenburg, Sweden

Keywords: SIMS, GCIB, Imaging, Depth-profiling, E. coli

Introduction

Antibiotic resistance causes 25 000 deaths every year. The rate of the spread and the threat that antibiotic resistance pose on a global scale makes this a field of great interest for research. Transfer of plasmids (mobile DNA) between bacterial cells through conjugation is one of the major contributors to the spread of antibiotic resistance. A larger understanding of the mechanisms behind this process and how to possibly hinder it could help to quench the development observed around the world. In this study, the changes in lipid composition in the cell membrane of three mutated E. coli strains that have previously been found unable to transfer plasmids were analysed using a J105-3D Chemical Imager.

Methods

The J105 (Ionoptika Ltd, UK) is an unconventional time-of-flight secondary ion mass spectrometry (ToF-SIMS) instrument specifically developed for analysis of organic materials [1]. A 40 keV CO₂ gas cluster ion beam (GCIB) is coupled to the J105 enabling imaging as well as depth profiling of native, underivatized biological samples and detection of intact lipids and other higher mass species [2]. Using the J105 coupled to a GCIB has allowed analysis of the outer membrane of mutated E. coli and how the lipid composition has changed in the membrane due to the mutations. Maximum autocorrelation factor (MAF) analysis was used as a tool to identify species of interest by mapping the largest differences within the data set.

Results

A previous study found that the stress response in E. coli exposed to carbon starvation were changes in the lipid composition by elongation, cyclopropanation, and increased cardiolipin formation [3]. In our study, three different mutated E. coli strains unable to transfer plasmids were analysed and compared to wild type E. coli. Here we have analysed mutants that have been shown to have impaired ability to directly transfer DNA (an important process in antibiotic resistance development). The sub-set of mutants selected for SIMS analysis possess mutations that are expected to influence lipid and protein composition in the cellular envelope of the bacteria and are thus ideally targets for surface sensitive MS analysis. The preliminary results from this study show extensive lipid changes in the mutated bacterial cell membranes, some of which are closely related to those previously observed in carbon starved E. coli cultures. Changes in elongation and cyclopropanation of unsaturated fatty acids were observed along with differences in the production of cardiolipins. Such effects can influence membrane fluidity and protein incorporation and also the functional properties of pili used to initiate conjugation. Lipid structure can also influence the fusion pore formation required for plasmid transfer during conjugation.

Conclusions

The ToF-SIMS GCIB/J105 system has here proven to be a good method for analysis of the surface and the sub-surface of the cellular envelope in mutated E. coli. MAF analysis was able to identify a large number of species which showed changes in the cell membrane of mutated strains relative to wild type E. coli. Some of these species have previously been identified as related to stress response in E. coli exposed to carbon starvation.

Novel Aspect

Linking lipid composition changes in the cell membrane in mutated E. coli to the bacteria's inability to transfer plasmids.

References

1. Fletcher, J. S., Rabbani, S., Henderson, A., Blenkinsopp, P., Thompson, S. P., Lockyer, N. P., Vickerman, J. C., A new dynamic in mass spectral imaging of single biological cells. *Analytical Chemistry*, 80, 9058-64 (2008)
2. Angerer, T. B., Blenkinsopp, P., and Fletcher, J. S., High Energy Gas Cluster Ions For Organic and Biological Analysis By Time-of-Flight Secondary Ion Mass Spectrometry. *International Journal of Mass Spectrometry*, 377, 591-598 (2015)
3. Wehrli, P. M., Angerer, T. B., Farewell, A., Fletcher, J. S., and Gottfries, J., Investigating the Role of the Stringent Response in Lipid Modifications during the Stationary Phase in *E. coli* by Direct Analysis with Time-of-Flight-Secondary Ion Mass Spectrometry *Analytical Chemistry*, 88, 8680-8688 (2016)

WP-160 / CHEMICAL IMAGING OF AGGRESSIVE BASAL CELL CARCINOMA USING TIME-OF-FLIGHT SECONDARY ION MASS SPECTROMETRY

Kelly Dimovska Nilsson⁽¹⁾ - *Marwa Munem*⁽¹⁾ - *Oscar Zaar*⁽²⁾ - *Noora Neittaanmäki*⁽²⁾ - *John S. Fletcher*⁽¹⁾ - *John Paoli*⁽²⁾

Department of Chemistry and Molecular Biology, University of Gothenburg, Gothenburg, Sweden (1) - Department of Dermatology and Venereology, University of Gothenburg, Gothenburg, Sweden (2)

Keywords: Imaging, SIMS, skin, cancer, GCIB

Introduction

ToF-SIMS has been successfully used for analysing lipids behaviour in biological samples like breast cancer tissue. Basal cell carcinoma (BCC) is one of the most increasing cancers worldwide and it is the most common malignancy in white people. Although the mortality is low as BCC rarely metastasises, this malignancy causes considerable morbidity and places a huge burden on healthcare services worldwide. Furthermore, people with BCC are at high risk of developing other malignancies.

Methods

Samples were collected from patients with basal cell carcinoma, by Mohs surgery. The tissue was cryo-sectioned for ToF-SIMS analysis and H&E staining of consecutive tissue slices was performed. ToF-SIMS was performed using an Ionoptika J105 instrument using a 40 keV (CO₂)6000+ ion beam.

Results

A set of basal cell carcinoma samples, removed by Mohs micrographic surgery and pathologically identified as having an aggressive subtype, have been analysed using time-of-flight secondary ion mass spectrometry (SIMS). The SIMS analysis employed a gas cluster ion beam (GCIB) to increase the sensitivity of the technique for the detection of intact lipid species. The GCIB also allowed these intact molecular signals to be maintained while surface contamination and delocalized chemicals were removed from the upper tissue surface. Distinct mass spectral signals were detected from different regions of the tissue (epidermis, dermis, hair follicles, sebaceous glands, scar tissue, and cancerous tissue) allowing mass spectral pathology to be performed. The cancerous regions of the tissue showed a particular increase in sphingomyelin signals that were detected in both positive and negative ion mode along with increased specific phosphatidylserine and phosphatidylinositol signals observed in negative ion mode. Samples containing mixed more and less aggressive tumor regions showed increased phosphatidylcholine lipid content in the less aggressive areas similar to a punch biopsy sample of a nonaggressive nodular lesion.

Conclusions

ToF-SIMS imaging with a GCIB distinguished cancerous and non-cancerous tissue from multiple patients along with different grades of high and low aggressive cancer.

Novel Aspect

ToF-SIMS analysis using 40 keV GCIB was used to analyse basal cell carcinoma samples, removed by Mohs micrographic surgery.

WP-161 / AN UNTARGETED APPROACH TO ELUCIDATE THE METABOLISM AND TISSUE DISTRIBUTION OF PHOSPHOROTHIOATE LINKED OLIGONUCLEOTIDES BY HIGH RESOLUTION MASS SPECTROMETRY

Bernd Steinhuber - Christophe Husser - Andreas Brink

F. Hoffmann-La Roche AG, PS, Basel, Switzerland

Keywords: MALDI-FTICRMS imaging, oligonucleotides, tissue drug disposition, target site ADME, metabolite identification

Introduction

Oligonucleotides linked with phosphorothioate bonds are a class of therapeutics under investigation in various pharmaceutical companies. Investigating the metabolism and disposition to target or non-target tissue by mass spectrometry could be a major advantage for a successful development. We developed two complementary analytical methods. To identify oligonucleotide metabolites an untargeted highly sensitive qualitative method using capillary flow liquid chromatography coupled to high resolution mass spectrometry was used, and for in situ imaging of tissue distribution a complementary method based on MALDI-FTICR MS. Both methods utilize a diagnostic fragment derived from the phosphorothioate backbone. This combined approach was successfully applied in a rat study with oligonucleotides of different sequence with and without conjugation to *N*-acetylgalactosamines (GalNAc).

Methods

Metabolite identification was achieved by capillary flow liquid chromatography with column switching coupled to a Fusion Orbitrap mass spectrometer. Data-dependent scans were performed in parallel to an untargeted MS² experiment (AIF) and a high resolution full scan. In situ tissue imaging was performed by MALDI-MS Imaging using a 7 Tesla Solarix XR FTICRMS instrument (Bruker). Matrix (α -Cyano-4-hydroxycinnamic acid) was applied using an iMatrixSpray device. Negative polarity and 300 laser shots were applied (2000 Hz) at 50 μ m lateral resolution. The quadrupole was used in full transferring mode followed by CID (60 eV). Full scan MS was acquired in the ICR cell. Tissue distribution of oligonucleotides was also evaluated using Immunohistochemistry.

Results

The detection and identification of oligonucleotides and their metabolites was achieved on a Fusion Orbitrap by combining full scan MS with two parallel MS² experiments, one data-dependent scan and an untargeted MS² experiment (AIF) applying high collision energy. In the AIF scan, a diagnostic fragment originating from the phosphorothioate backbone (O_2PS^- : m/z 94.936) was formed efficiently upon collisional activation. Based on this fragment an accurate determination of oligonucleotide metabolites was achieved, independent of their sequence in an untargeted but highly selective manner. The method was effectively applied to investigate uptake and metabolism of oligonucleotides in rat liver and rat kidney homogenates down to the subnanomolar range. For MALDI-FTICR MS Imaging, oligonucleotides were detected indirectly and independently of their sequence using the quadrupole in an unselective transferring mode followed by strong collisional activation. We could identify in the ICR cell the diagnostic fragment (O_2PS^- : m/z 94.936). Both methods were applied in a mechanistic multiple-dose rat study. Oligonucleotides with different sequence were used both as unconjugated compounds or conjugated with *N*-acetyl galactosamine which is a high-affinity ligand for the asialoglycoprotein receptor on the hepatocyte cell surface. Liver and kidney were investigated for metabolites and were evaluated in situ by immunohistochemistry and MALDI-FTICRMS Imaging. In the liver, unconjugated oligonucleotides were

demonstrated mainly in Kupffer cells, whereas GalNAc-conjugated oligonucleotides accumulated in both Kupffer cells and hepatocytes. In the kidney, all oligonucleotides tested accumulated in the cortex primarily in proximal tubules. Interference with nucleotides of endogenous RNAs or DNAs was not observed. Phosphorothioate instead of phosphate linkage is a common modification for therapeutic oligonucleotides. Therefore, this approach might be generally applicable to investigate tissue distribution of this evolving therapeutic modality.

Conclusions

It could be shown that the distribution and metabolism of oligonucleotides linked with phosphorothioate bonds can be analyzed in an untargeted manner by mass spectrometry. The metabolism of oligonucleotides was elucidated by LC-MS/MS whereas in situ imaging of tissue distribution was achieved on a MALDI-FTICR MS instrument. Both approaches harness the phosphorothioate backbone of these oligonucleotides as a diagnostic fragment (O_2PS^- : m/z 94.936) as this is not present in endogenous nucleic acids. Both methods proved their applicability in a mechanistic multiple-dose rat study for which several oligonucleotides of different sequences were dosed.

Novel Aspect

Untargeted approach to investigate metabolism and tissue distribution of therapeutic phosphorothioate linked oligonucleotides by high resolution LC-MS and MALDI-FTICRMS Imaging.

References

1. Husser C., Brink A., Zell M., Müller M.B., Koller E., Schadt S., *Anal. Chem.*, 89, 6821–6826 (2017)

WP-162 / A CASE STUDY FOR DEMONSTRATING THE USEFULNESS OF IMAGING MASS SPECTROMETRY IN TOXICITY ASSESSMENT USING AN AMIODARONE-INDUCED PHOSPHOLIPIDOSIS MODEL

Koji Tanaka⁽¹⁾ - ***Akane Kashimura***⁽²⁾ - ***Hiroko Sato***⁽²⁾

Mitsubishi Tanabe Pharma Corporation, DMPK Research Laboratories, Toda, Japan (1) - Mitsubishi Tanabe Pharma Corporation, Safety Research Laboratories, Toda, Japan (2)

Keywords: imaging mass spectrometry, toxicity assessment, histological examination, amiodarone

Introduction

Investigating the local distribution of drug and biomarker in target organs provides useful information, which accelerates drug discovery and development. Imaging mass spectrometry (IMS) is a powerful technique to detect the localization of compounds of interest in biological tissue without any labeling. In the present report, we demonstrate the utility of IMS to investigate drug and biomarker distribution in an amiodarone-induced phospholipidosis model.

Methods

Rats were orally administered amiodarone (AMD) for 7 days at 150 mg/kg/day. After the final dose, lungs, spleens and mesenteric lymph nodes were collected. Prepared cryosections (5µm/section) were coated with α -cyano-4-hydroxycinnamic acid using iMLayer (Shimadzu) and analyzed using the IMS system (Q-Exactive (Thermo Scientific) equipped with MALDI ion source AP-SMALDI 10 (TransMIT)). Mass images were generated using Million ver. 3.2 software (TransMIT).

Results

AMD, its metabolites (*N*-desethyl amiodarone and mono-oxidized amiodarone), and representative phospholipids, which are considered biomarkers of phospholipidosis, were detected in lungs, spleens and mesenteric lymph nodes in a single IMS analysis. Detection intensity of all the compounds was relatively

higher in foamy macrophage infiltration areas than in non-infiltration areas. For the spleen, the distribution pattern of representative phospholipids markedly differed between AMD-administered rat and vehicle-administered rat.

Conclusions

These results indicated that AMD and its metabolites were stored in macrophages, and this is quite consistent with the mechanism of amiodarone-induced phospholipidosis. They also indicate that IMS is a useful technique for evaluating the distribution of drugs and biological components in the elucidation of toxicity mechanisms.

Novel Aspect

IMS analysis revealed the distribution pattern of AMD, its metabolites, and representative phospholipids in an amiodarone-induced phospholipidosis model.

References

1. Akane K., Kouji T., Hiroko S., Hidefumi K., Masaharu T., *Journal of Toxicologic Pathology*, 31, advance online publication (2018).

WPS-S07 - Gas phase ion chemistry & thermochemistry

Chair: Pan Yanjiang

Wednesday

WP-163 / A QUANTUM CHEMICAL TOOL FOR THE PREDICTION OF MASS SPECTRA AND THE IDENTIFICATION OF FRAGMENTATION PATHWAYS

Julie Cautereels - Frank Blockhuys

University of Antwerp/Departement of Chemistry, Antwerp

Keywords: Quantum Chemical Mass Spectrometry (QCMS2), Density Functional Theory (DFT), inter-side-chain (ISC) interactions, Tandem Mass Spectrometry, Tripeptides

Introduction

A new computational tool for the prediction of mass spectra based on quantum chemical calculations is being developed, called Quantum Chemical Mass Spectrometry for Materials Science (QCMS2). It was benchmarked by predicting the EI fragmentation pathways of a number of organics with a variety of functionalities. The main features in the mass spectra were correctly reproduced and new fragmentation routes proposed by QCMS2 were confirmed by MS/MS experiments [1].

Methods

QCMS2 is based on DFT/B3LYP/6-311+G* calculations and selects fragmentation pathways based on the relative values of bond order, reaction energies (for bond cleavages) and transition state energies (for rearrangements). Known peptide-specific fragmentation mechanisms such as the Mobile Proton Model (MPM) for His, confirmed by quantum chemical calculations, have been implemented into the method and the code omits the breaking of irrelevant bonds [2].

Results

QCMS2 is being applied to predict the fragmentation patterns of tripeptides in ESI CID MS/MS focusing on the influence of inter-side-chain (ISC) interactions on the initial protonation and subsequent fragmentation.

The fragmentations of a number of non-cyclic tripeptides consisting of His as the central amino acid (AA) and AAs which allow the formation of strong hydrogen bonds as the peripheral ones, have been studied. Observations such as (1) the fact that the (dis)appearance of ISC interactions in the Lys-His-X (X = Asn, Asp, Gln, Glu and Ser) and Trp-His-X combinations determines the most stable protonation site of the tripeptide, and (2) the fact that in Arg-His-Ser and Asn-His-Ser ISC interactions lead to specific fragmentations, illustrate the importance of strong interactions between AA side chains. The results obtained from QCMS2 are compared to those obtained from conventional methods to demonstrate how many of the peaks not predicted by these latter methods are due to fragmentations involving the side chains and are predicted by QCMS2.

Conclusions

QCMS2 correctly reproduces the main features in the mass spectra of tripeptides, *i.e.*, the traditional backbone cleavages resulting in α -, β 2 and γ 1-ions. Furthermore, due to its taking the fragmentation of the side chain and ISC interactions into account, it also predicts fragments resulting from these and, as a result, convincingly outperforms common tools such as PEAKS, PeptideArt and Prospector.

Novel Aspect

QCMS2 assigns up to twice the number of signals in peptide mass spectra assigned by common mass spectral assignment software, and generates detailed fragmentation routes.

References

1. Cautereels J., Claeys M., Geldof D., Blockhuys F., *J. Mass Spectrom.*, 51, 602 (2016).
2. Cautereels J., Blockhuys F., *J. Am. Soc. Mass Spectrom.*, 28, 1227 (2017).

WP-164 / ISOMERIZATION AND DISSOCIATION OF ALPHA-CENTERED AND PI-CENTERED GLYCYLGLYCYLTRYPTOPHAN RADICAL CATIONS: STRUCTURAL, ENERGETIC, MECHANISTIC, AND SPECTROSCOPIC INVESTIGATIONS

Andy C. K. Siu⁽¹⁾ - **Jonathan Martens**⁽²⁾ - **Andre Sze**⁽³⁾ - **Jos Oomens**⁽²⁾ - **Ivan K. Chu**⁽³⁾

City University of Hong Kong, Department of Chemistry, Hong Kong (1) - Institute for Molecules and Materials, FELIX Laboratory, Nijmegen (2) - The University of Hong Kong, Department of Chemistry, Hong Kong (3)

Keywords: carbon-centered radicals, glycyglycyltryptophan, gas-phase

Dissociation of ligated copper(II)–peptide complexes, $[\text{CuII}(\text{L})\text{M}]^{\bullet 2+}$ (L, ligand; M, peptide), generates peptide radical cations ($\text{M}^{\bullet +}$) through single electron transfer. [1] Diversity is introduced into the ion chemistry of a peptide radical cation through the possibility of the radical migrating in a number of different ways to corresponding prerequisite radical sites prior to fragmentation; some of these competitive reactions are charge- and radical carrier-specific, and can be perturbed upon varying the nature or location of the radical/charge site, the basicity of the residue at that site, or the peptide sequence, all of which can substantially impact the isomerization barriers for radical or charge mobility and the subsequent reactions. [2–5] The major impetus of this study was to use a simple prototypical system, $[\text{CuII}(\text{L})\text{YGW}]^{\bullet 2+}$, to track the location of the radical and charge sites, starting from the copper(II)–peptide complex to the predetermined radical sites of the peptide radical cations after a complex permutation of intramolecular hydrogen atom migrations—potentially allowing us to determine whether or not the individual isomerization, rearrangement, or fragmentation pathways would predominate. A wide variety of site-specific radical peptides can be produced efficiently with judicious choice of the metal–ligand–peptide system; [1–2] we used isomeric π -centered $[\text{GGW}\pi]^{\bullet +}$ (formed directly through dissociation of $[\text{CuII}(\text{amine})\text{GGW}]^{\bullet 2+}$ complexes) and α -centered $[\text{G}\bullet\text{GW}]^+$ (generated through side-chain cleavage of the tyrosyl residue of the

YGW^{•+} radical cation) species for our initial infrared multiple photon dissociation (IRMPD) investigations, in conjunction with density functional theory calculations. The IRMPD spectrum of G•GW⁺ displays a strikingly different set of product ions than those obtained from [GGWπ]^{•+}. It appears that the two isomeric radical cationic GGW structures did not interconvert on the time scale of our IRMPD experiments. Such isomeric π- and σ-centered radical cationic peptides comprise a well-defined prototypical system for dissecting the fundamental parameters that govern the isomerization and dissociation of peptide radical cations, which undergo various rearrangement processes, including hydrogen atom, proton, and electron transfer, in the absence of solvation.

References

1. Chu I. K., Rodriguez C. F., Lau T. C., Hopkinson A. C., Siu K. W. M., *J. Phys. Chem. B*, 104, 3393–3397 (2000)
2. Barlow C. K., McFadyen W. D., O'Hair R. A. J., *J. Am. Chem. Soc.*, 127, 6109–6115 (2005)
3. Hopkinson A. C., Siu K. W. M., *Principles of Mass Spectrometry Applied to Biomolecules*, Laskin J., Lifshitz C., Eds., John Wiley and Sons: New York (2006)
4. Chu I. K., Laskin J., *Eur. J. Mass Spectrom.*, 17, 543–556 (2011)
5. Turecek F., Julian R., *Chem. Rev.*, 113, 6691 (2013)
6. Hopkinson A. C., *Mass Spectrom. Rev.*, 28, 655–671 (2009)
7. Mu X., Song T., Siu C. K., Chu I. K., *Chem. Record*, 17, 1–26 (2017)
8. Chu I. K., Hopkinson, A. C., Siu K. W. M., *J. Am. Chem. Soc.* 2008, 130, 7862–7872
9. Laskin J., Futrell J. H., Chu I. K., *J. Am. Chem. Soc.*, 129, 9598–9599 (2007)
10. Ng C. M., Song T., Laskin J., Chu I. K., *J. Phys. Chem. B*, 114, 2270–2280 (2010).
11. Easton C. J. *Chem. Rev.*, 97, 53–82 (1997)
12. Bordwell F. G., Lynch T. Y. *J. Am. Chem. Soc.*, 111, 7558–7562 (1989)
13. Viehe H. G., Janousek Z., Merenyi R., Stella L., *Acc. Chem. Res.*, 18, 148–154 (1985).
14. Simon, S., Sodupe M., Bertran J., *J. Phys. Chem. A*, 106, 5697–5702 (2002).
15. Simon S., Gil A., Sodupe M., Bertran J., *J. Mol. Struct., Theochem.*, 727, 191–197 (2005)
16. Polce M. J., Wesdemiotis C., *J. Am. Soc. Mass Spectrom.*, 10, 1241–1247 (1999).

WP-165 / NITROGEN OXIDES BINDING TO COBALAMIN: ION CHEMISTRY BY FT-ICR MASS SPECTROMETRY

Maria Elisa Crestoni - Alessandro Maccelli - Simonetta Fornarini

Sapienza Università di Roma, Dipartimento di Chimica e Tecnologie del Farmaco, Roma

Keywords: ESI-FT-ICR, Cobalamin, nitric oxide, ion-molecule reactions

Introduction

Cobalamines (Cbls) are important cofactors in various enzymatic processes. Many studies indicate that the reduced form of cobalamin, Cbl(II), reacts reversibly with NO to form a nitrosylcobalamin adduct at a high rate. This evidence is supported by the NO-induced inactivation of Cbl-dependent enzymes and the inhibition of NO actions by Cbl supplementation. [1-2] Here, the intrinsic features for NO_x (x=1,2) binding to bare Cbls ions are assessed by ESI FT-ICR MS.

Methods

Experiments were carried out with a Bruker BioApex 4.7 FT-ICR mass spectrometer coupled with an ESI source. The reduced Cbl(II) and oxidized Cbl(III) forms were formed by ESI of a 10 μM aqueous hydroxycobalamin solution at physiological pH, mass-selected and exposed to NO_x (x = 1, 2) in the ICR cell at stationary pressures of 1.0-10 × 10⁻⁸ mbar. Pseudo-first order-rate constants divided by the neutral concentration afford second-order rate constants. [3]

Results

In order to improve the current understanding of NO_x (x = 1, 2) scavenging by Cbls, the kinetics and mechanistic behavior of several Cbls species, either in the reduced or oxidized and protonated/sodiated or deprotonated forms, have been explored in the gas phase, lacking any perturbation by solvent or counterions. Some relevant/preliminary evidences may be outlined: i) both sampled neutrals undergo an addition process; ii) for NO₂ binding, an overall higher yield is observed with either Cbl(II) or Cbl(III), by a factor of ca. 3 and 10, respectively. The efficiency values span from 4.0 to 18% for Cbl(III) and doubly deprotonated Cbl(III), the only form which undergoes a double addition, whereas (de)protonated and monosodiated Cbl(II) show similar values (~ 8%); iii) the reactivity with NO shows a clear bias for Cbl(II) species; iv) Cbl(III) species undergo also hydrogen atom abstraction by NO_x (x = 1, 2).

Conclusions

Gas phase ion-molecule reactions provide mechanistic information relevant to biological functions void of environmental effects. Here, the binding properties of nitrogen oxides to various forms of cobalamin have been addressed, involving either deprotonated "base-on" or protonated "base-off" species. This study confirms preferential NO ligation to Cbl(II) and outlines the influence of the coordination environment of the cobalt ion.

Novel Aspect

Direct determination of reaction patterns and thermal rate constants are provided for NO_x (x = 1, 2) binding to naked Cbl(II) and Cbl(III) ions in a solvent free environment.

References

1. Dereven'kov I. A., Salnikov D. S., Silaghi-Dumitrescu R., Makarov S. V., Koifman O. I., *Coord. Chem. Rev.* 309, 68-83 (2016).
2. Selçuki C., van Eldik R., Clarck T., *InorgChem*, 43, 2828-2833 (2004).
3. Cantú Reinhard F. G., Sainna M. A., Upadhayay P., Balan G. A., Kumar D., Fornarini S., Crestoni M. E., de Visser S. P., *Chem. Eur. J.* 22, 18608-18619 (2016).

WP-166 / TANDEM MASS SPECTROMETRY AS A TOOL TO EVALUATE LIGAND ELECTRONIC EFFECT ON AU(I) COMPLEXES

David Gatineau⁽¹⁾ - **Denis Lesage**⁽²⁾ - **Héloïse Dossmann**⁽²⁾ - **Hervé Clavier**⁽³⁾ - **Anne Milet**⁽¹⁾ - **Yves Gimbert**⁽¹⁾

DCM UMR 5250, CNRS-Université Grenoble Alpes, Grenoble (1) - IPCM, Sorbonne Université, CNRS, Paris (2) - ISM2, Aix Marseille Université CNRS, Centrale Marseille, Marseille (3)

Keywords: Bond dissociation energy, RRKM modelling, CID, tandem mass spectrometry, gold chemistry

Introduction

Homogenous gold catalysis is a powerful approach to access high diversity and complexity of molecules. [1] Besides the metal, ligands play a major role to control and orientate the reactivity. [2] Many theoretical works have aimed and succeeded to interpret the particular behavior of ligands in cationic gold complex. [3] In contrast, experimental approaches to evaluate these effects and validate these theoretical models are scarce and would be highly desirable.

Methods

A collection of gold carbonyl complexes was studied using a modified triple-quadrupole equipped with an ESI source. L-Au-CO⁺ ions were directly formed in the source and submitted to Collision Induced Dissociation (CID) experiments. From survival yield (SY) breakdown curves, critical energy values were obtained using a kinetic model based on transition state theory (RRKM theory). results were analyzed using

DFT calculations and energy decomposition analysis (EDA).

Results

Au(I) complexes were successfully formed in the ESI source for 16 L ligands: expected, the greater the capacity of L to enrich gold, the higher the dissociation energy of Au-CO. The SY for this dissociation were recorded for different collision energies during CID experiments. Due to the variety of size and masses of the complexes, a semi-empirical calibration had to be established in order to get an evaluation of the mean internal energies of the complexes prior to decomposition. From there, critical energies were obtained using kinetic modelling. The variety of electronic effects of the studied ligands is reflected by these values: a range of 10.5 kcal/mol is observed between the weakest and the strongest Au-CO bond, as confirmed by DFT calculations. To gain insights into the Au-CO bond nature, an energy decomposition analysis was also performed. The contributions of the orbital and electrostatic interactions and the Pauli repulsion to the bonding were considered. results are rationalized by considering the influence of each of these terms.

Conclusions

Bond dissociation energies of 16 gold-carbonyl complexes were obtained using mass spectrometry and kinetic modelling. Depending on the nature of the ligand, a more or less large difference in the BDEs could be observed. These data were assessed as new probes for evaluating ligand effects in gold catalysis. With the help of energy decomposition analysis, a rationalization of the various electronic effects observed for the studied ligands could be achieved.

Novel Aspect

A new experimental approach, based on mass spectrometry experiments and kinetic modelling, is proposed to evaluate the ligands effects in gold(I) complexes.

References

1. Dorel R., Echavarren A. M., Chem. Rev. 115, 9028. (2015).
2. Gorin D. J., Sherry B. D., Toste F. D., Chem. Rev., 108, 3351 (2008).
3. Bistoni G., Rampino S., Scafuri N., Ciancaleoni G., Zuccaccia D., Belpassi L., Tarantelli F., Chem. Sci. 7, 1174 (2016).

WP-167 / ENERGY-RESOLVED COLLISION-INDUCED DISSOCIATION STUDY OF PROTON-BOUND HETRODIMERS OF GUANINE AND CYTOSINE

Sang Yun Han⁽¹⁾ - Jeong Ju Park⁽¹⁾ - Choong Sik Lee⁽²⁾

Gachon University, Department of Nanochemistry, Seongnam (1) - Ministry of National Defense, Scientific Investigation Lab., Seoul (2)

Keywords: Proton-bound base pairs; Guanine; Cytosine; ER-CID; DFT

Introduction

To study the effect of protonation on DNA base pairing, we investigated energy-resolved collision-induced dissociation (ER-CID) of proton-bound heterodimers of guanine (G), guanosine (Gn), cytosine (C), and 1-methylcytosine (1-MeC), and cytidine (Cn), which are likely to be formed by Hoogsteen pairing including an ionic hydrogen bond. The structures, energies, and proton transfer reactions of dimeric conformers were also investigated theoretically.

Methods

The protonated dimers of $[G:C:H]^+$, $[G:1-MeC:H]^+$, and $[Gn:Cn:H]^+$ were produced by electrospray ionization in the acidified conditions, which were interrogated by ER-CID experiments using LTQ XL (Thermo Fisher) and Xevo TQ (Waters) for multiple and single collision conditions, respectively. For quantum chemical calculations, density functional theory (DFT) at the B3LYP/6-311+G(2d,p) level using Gaussian 09 program

suite was utilized.

Results

Recently, an extensive theoretical study[1] and an IR multiple photon dissociation (IRMPD) spectroscopy[2] suggested that the Hoogsteen basepair, *i.e.* CH⁺•••G, be the dominant conformation for [G:C:H]⁺ in the gas phase. Based on the previous investigations, we modelled the proton-bound heterodimers of [G:C:H]⁺, [G:1-MeC:H]⁺, and [Gn:Cn:H]⁺ as Hoogsteen base pairs using DFT calculations. ER-CID behaviours of the protonated heterodimers were explored both in multiple and single collision conditions. It was found that dissociation of [G:C:H]⁺ led to preferential formation of C:H⁺, while simple consideration of the kinetic method [3,4] predicts preferential protonation on more proton-loving moiety, G. However, CID of the other basepairs of [G:1-MeC:H]⁺ and [Gn:Cn:H]⁺ occurred with favored protonation on more proton-loving moieties, yielding 1-MeC:H⁺ and Gn:H⁺ as the dominant fragments, respectively.

Conclusions

DFT and CID study suggested Hoogsteen structures for [G:C:H]⁺, [G:1-MeC:H]⁺, and [Gn:Cn:H]⁺. By varying substituents, alternation of energetics in base-pairing, dissociation, and intermolecular proton transfer reactions were introduced and examined. While the three heterodimers have similar energetic relations, CID of [G:C:H]⁺ showed a unique behavior, where certain dynamic reasons are likely to be deeply involved.

Novel Aspect

ER-CID of gas-phase Hoogsteen base pairs occurring by protonation of bases was examined for the first time, and an intriguing dissociation of [G:C:H]⁺ was revealed.

References

1. Jun J., Han S. Y., Theo. Chem. Acc., 136, 136 (2017).
2. Cruz-Ortiz A. F., Rossa M., Berthias F., Berdakin M., Maitre P., Pino G. A., J. Phys. Chem. Lett., 8, 5501 (2017).
3. Cooks R. G., Wong P. S. H., Acc. Chem. Res., 31, 379 (1998).
4. Seong Y., Han S. Y., Jo S. -C., Oh H. B., Mass Spectrom. Lett., 2, 69 (2011).

WP-168 / CORRELATION BETWEEN STRUCTURE AND STABILITY OF GAS PHASE OLIGONUCLEOTIDES PROBED BY TUNED COLLISION- AND VUV-MS/MS

Stefanie Ickert⁽¹⁾ - Jens Riedel⁽²⁾ - Sebastian Beck⁽³⁾ - Michael W. Linscheid⁽³⁾

Humboldt Universität zu Berlin, Federal Institute for Materials Research and Testing, Berlin (1) - Federal Institute for Materials Research and Testing, Berlin (2) - Humboldt Universität zu Berlin, Berlin (3)

Keywords: DNA, Locked Nucleic Acids, VUV dissociation/postionization, CID fragmentation, Ion Mobility

Introduction

Various MS/MS methods are commonly used to create comprehensive information for bioanalytical applications. However, only little data is available on gas phase and fragmentation behavior of important molecules like double stranded DNA and related clinically relevant analogs. Previously, CID was predominantly used to characterize these kind of samples, leaving other techniques underused. [1]

Methods

A LTQ Orbitrap XL was used to generate CID and HCD spectra of single and double strand DNA and DNA-LNA hybrids. The Synapt G2-S was used for collisional cross section assessment. As an additional MS/MS method a continuous vacuum ultraviolet lamp with transmission maxima at 125 and 160 nm was used. This approach provided efficient fragmentation along atypical, orthogonal pathways in both ionization modes

and sensitivity increase by postionization of neutral losses.

Results

DNA folding in the gas phase is strongly dependent on the charge level. With the use of IMS and CID data, for this unfolding a clear correlation between the energetic dissociation onset and the reorganization could be determined. [1,2] LNA, as a more rigid structure in comparison to DNA, indicated a complete different behavior concerning unfolding and fragmentation pathways. This was not only found for single strand nucleic acids, but also for the respective DNA-DNA and DNA-LNA duplexes. While DNA dissociates along energy specific pathways, LNA was found to requires higher energy for fragmentation and follow less specific dissociation channels.

Additionally, VUV fragmentation was tested for positive and negative ionization. [3] For multiply charged oligonucleotides and proteins a charge reduction following electron ejection followed by dissociation was found. Samples of other categories showed fragmentation into common and atypical fragments. Moreover, for some systems, a strong partition of radiation induced postionization was observed.

Conclusions

Double and single strands of DNA and LNA were selectively fragmented. While DNA showed clear pathways including reorganization, the rigidity of LNA caused multiple dissociation channels and different folding pattern. VUV as new MS/MS method produced charge reduction of multiple charged ions as well as new fragment classes. Also an increase of intensity was observed by postionization of fragment ions.

Novel Aspect

A comprehensive study of DNA and LNA concerning gas phase organization and dissociation in MS and IMS was achieved. VUV was introduced as a new MS/MS method.

References

1. S. Ickert, J. Riedel, S. Beck, M. W. Linscheid, Fragment Formation of Oligonucleotides in Collision Induced Mass Spectrometry Experiments, *Rapid Commun. Mass Spectrom.*, 32, 597 (2018)
2. S. Ickert, J. Hofmann, J. Riedel, S. Beck, K. Pagel, M. W. Linscheid, From Solution to Isolation: Non covalent Interactions Stabilize Compact Structures of DNA Oligomers in the Gas Phase, *Eur. J. Mass Spectrom.* 24, 225 (2018)
3. S. Vogel, J. Rackwitz, R. Schürmann, J. Prinz, A. Milosavljevic, M. Refregiers, A. Giuliani, I. Bald, Using DNA Origami Nanostructures to Determine Absolute Cross Sections for UV Photon-Induced DNA Strand Breakage, *J. Phys. Chem. Lett.*, 6, 4589 (2015)

WP-169 / MEASUREMENT OF PROTONATED WATER CLUSTERS USING HIGH KINETIC ENERGY ION MOBILITY SPECTROMETRY (HIKE-IMS)

Ansgar T. Kirk - Maria Allers - Stefan Zimmermann

Leibniz Universität Hannover, Department of Sensors and Measurement Technology, Hannover

Keywords: protonated water cluster, hydronium ion, cluster distribution, ion mobility spectrometry

Introduction

Protonated water clusters $(\text{H}_2\text{O})_n\text{H}^+$ play an important role during many ionization processes, including atmospheric pressure chemical ionization (APCI), and in particular proton transfer reaction (PTR). The cluster size n influences ionization by changes in proton affinity and ion mobility. This is a major issue when calculating concentrations in PTR-MS [1, 2]. Furthermore, transfer into vacuum changes the water cluster distributions, rendering them unknown.

Methods

High Kinetic Energy Ion Mobility Spectrometry (HiKE-IMS) ionizes in a reaction tube operated at reduced pressures (20 mbar) and high reduced field strengths (up to 130 Td) [3], similar to PTR-MS. However, instead of a mass spectrometer, an ion mobility spectrometer operated at the same pressure as the reaction tube is used for ion separation and detection. This allows analyzing the ions under the ionization conditions with anion mobility resolving power of $R_p=140$ [4].

Results

When monitoring the ion mobility spectrum while sweeping the reduced field strength in both the reaction and the drift region, a shift of the protonated water cluster peak position can be observed at the energies where a change in most prevalent cluster size is expected (30 Td, 60 Td, 90 Td). Significantly below each of these thresholds, a Gaussian peak was measured at the expected mobility of the respective larger water cluster. Near the threshold, the peak starts broadening and a changing ion distribution between the positions of the respective larger and smaller water cluster can be observed. Significantly above each threshold, again a Gaussian peak was observed, this time at the position of the smaller water cluster.

Conclusions

Using HiKE-IMS allows analyzing protonated water clusters under the conditions inside the ionization source. This way, both the prevalent cluster size and the ion mobility can be determined as for example needed for correct concentration calculation in PTR-MS.

Novel Aspect

For the first time, protonated water cluster distributions and ion mobilities have been directly measured under the conditions present during ionization.

References

1. Warneke C., van der Veen C., Luxembourg S., de Gouw J. A., Kok A., *Int. J. Mass Spec.*, 207, 167-182 (2001)
2. Keck L., Oeh U., Hoeschen C., *Int. J. Mass Spec.*, 264, 92-95 (2007)
3. Langejuergen J., Allers M., Oermann J., Kirk A. T., Zimmermann S., *Anal. Chem.*, 86, 11841-11846 (2014)
4. Kirk A. T., Grube D., Kobelt T., Wendt C., Zimmermann S., *Anal. Chem.*, DOI: 10.1021/acs.analchem.7b04586

Acknowledgements

Funded by the Deutsche Forschungsgemeinschaft (DFG, German Research Foundation) – ZI 1288/7-1.

WP-170 / CHALLENGES OF SOFT CHEMICAL IONIZATION MASS SPECTROMETRY: ION CHEMISTRY OF GLYOXAL

Michal Lacko⁽¹⁾ - ***Felix Piel***^(2,3) - ***Andreas Mauracher***⁽³⁾ - ***Antonín Knížek***⁽¹⁾ - ***Patrik Španěl***⁽¹⁾

J. Heyrovský Institute of Physical Chemistry of the CAS, Prague (1) - IONICON Analytik GmbH (2) - Institut für Ionenphysik und Angewandte Physik, Leopold-Franzens-Universität Innsbruck, Innsbruck (3)

Keywords: ion molecular kinetic, soft chemical ionization, glyoxal, SIFT-MS

Introduction

Glyoxal (m. w. 58 g/mol) is a highly reactive molecule associated with some basic biological processes [1], photo-catalysis of CO_2 into CH_4 on TiO_2 or photo-dissociation of anthropogenic and biogenic volatile organic compounds (VOC) [3]. Presence of glyoxal can be monitored using soft chemical ionization mass spectrometry (SCI-MS). Here we discuss the ion chemistry of glyoxal interacting with H_3O^+ , NO^+ and $\text{O}_2^{+\bullet}$ reagent ions used in SCI-MS and possibility of detection.

Methods

Selected ion flow tube mass spectrometry, SIFT-MS[4], was used to investigate ion-molecule reactions of H_3O^+ , NO^+ and O_2^{*+} reagent ions with glyoxal at variable humidity. The ion chemistry under thermal conditions is theoretically described using B3LYP/6-311(p,d) calculations and using numerical modelling of the ion-molecule reaction kinetics.

Results

NO^+ reactions lead to formation of $\text{C}_2\text{H}_2\text{O}_2(\text{NO})^+$ and $\text{C}_2\text{H}_2\text{O}_2^{*+}$ ions, in the humid environment the association reaction is dominant. Charge transfer was observed for O_2^{*+} reactions as well. Additional dissociation observed previously for O_2^{*+} reagent ion was not confirmed [5-6]. The H_3O^+ reaction proceeds at dry conditions via proton transfer reaction forming $\text{C}_2\text{H}_3\text{O}_2^+$. With increasing concentration of water molecules, formation of m/z 31 (CH_2OH^+) is becoming dominant. Corresponding hydrated ions were observed as well. Using quantum chemistry calculations, we estimated $\text{PA}(\text{C}_2\text{H}_2\text{O}_2) = 165.1$ kcal/mol and characterised the energetics of reaction pathways forming CH_2OH^+ . Reaction kinetics modelling indicates formation of a stable intermediate $\text{C}_2\text{H}_3\text{O}_2^+ \cdot \text{H}_2\text{O}$ which dissociates to protonated formaldehyde CH_2OH^+ after interaction with water.

Conclusions

According to theory, protonated formaldehyde is formed from $\text{C}_2\text{H}_3\text{O}_2^+ \cdot \text{H}_2\text{O}$ by reaction with water resulting into $\text{CH}_2\text{OH}^+ + \text{HCOOH} + \text{H}_2\text{O}$ products or by a formation of excited $\text{C}_2\text{H}_3\text{O}_2^+ \cdot \text{H}_2\text{O}^*$ after two body association which dissociates to CH_2OH^+ directly. Finally, protonated formaldehyde hydrate is not the terminating ion at high humidities. Thus, we expect a catalytic effect of $\text{CH}_2\text{OH}^+ \cdot \text{H}_2\text{O}$ forming a proton bound water dimer as $\text{CH}_2\text{OH}^+ \cdot \text{H}_2\text{O} + \text{H}_2\text{O} \rightarrow \text{CH}_2\text{O} + \text{H}_3\text{O}^+ \cdot \text{H}_2\text{O}$.

Novel Aspect

Understanding of ion chemistry of glyoxal is important for its analysis in gas phase using SCI-MS techniques and its distinguishing from isobaric compounds.

References

1. Manini P., d'Ischia M., As Markers of Redox Imbalance in Fanconi Anemia and in Other Oxidative Stress-Related Disorders. Fanconi Anemia And Oxidative Stress, Nova Science Publishers, New York 2015, p. 69.
2. Shkrob I. A., Marin T. W., He H. and Zapol P., J. Phys. Chem. C 116, 9450-9460 (2012).
3. Volkamer R., Platt U., Wirtz K., J. Phys. Chem. A 105, 7865 (2001).
4. Smith D., Španěl P., Holland T. A et al., Rapid Commun. Mass Spectrom. 13, 724 (1999).
5. Michel E., Schoon N. et al., Int. J. Mass Spectrom. 244, 50-59 (2005).
6. Guimbaud C., Catoire V. et al., Int. J. Mass Spectrom. 263, 276-288 (2007).

WP-171 / WALK-ON-SPHERE REARRANGEMENT AND RETRO-BINGEL REACTION OF GAS-PHASE FULLERENE MALONATE IONS

Johannes Oswald - Martin B. Minameyer - Katerina L. Maxouti - Andreas Hirsch - Thomas Drewello

Friedrich-Alexander-Universität Erlangen-Nürnberg, Department Chemie und Pharmazie, Erlangen

Keywords: APCI, organic fullerene derivatives, retro-Bingel, fragmentation mechanism

Introduction

C_{60}LX ($\text{L}=\text{C}(\text{COOEt})_2$; $x=4-6$) are investigated by APCI with a QIT- and a QTOF-MS to study their chemical stability as radical ions in the gas phase. For uncharged (neutral) fullerene adducts, the positions of the cyclopropane rings on the C_{60} surface are highly stable. In contrast, under electrochemically reducing conditions in solution an isomerization of the molecules via a walk-on-sphere rearrangement is observed [1].

Methods

APCI experiments were performed with a QIT- and a QTOF-MS. The samples were synthesized by a published procedure [2]. C₆₀L₆ and C₆₀L₅ are isomerically pure (determined by ¹³C NMR). C₆₀L₆ has an octahedral addition pattern with TH symmetry. C₆₀L₅ has a quadratic pyramidal addition pattern with C_{2v} symmetry. C₆₀L₄ consists of a mix of isomers with varying symmetries.

Results

For all examined fullerene adducts, radical anion (C₆₀LX^{•-}) and radical cation (C₆₀LX^{•+}) formation is observed applying APCI. Fragmentation experiments with a QIT-MS show the repeated loss of L₂ as the main fragmentation pathway for C₆₀LX^{•-} and C₆₀LX^{•+} (x=4-6). [3] Energy-dependent fragmentation experiments using a QTOF-MS show a decreasing energy-demand for the fragmentation when the number of cyclopropane moieties increases. Furthermore, dissociation of radical anions is less energy demanding than of radical cations, this observation is in agreement with experiments by Kessinger et al. which show isomerization of C₆₀L₂ in solution under electrochemically reducing conditions by a walk-on-sphere rearrangement [1].

Conclusions

Both C₆₀LX^{•-} and C₆₀LX^{•+} (x=4-6) show a very selective fragmentation with the loss of L₂. This fragmentation can only occur when two diethyl ligands "walk" on the surface towards each other. The walk-on-sphere rearrangement is followed by the formation of a malonate dimer (L₂) which is then released from the fullerene.

Novel Aspect

The charge-induced isomerization of C₆₀LX (L=C(COOEt)₂; x=4-6) has been demonstrated unequivocally by observation of the loss of a malonate dimer (L₂).

References

1. R. Kessinger, M. Gómez-López, C. Boudon, J. Gisselbrecht, M. Gross, L. Echegoyen, F. Diederich, J. Am. Chem. Soc. 120, 8545 (1998)
2. A. Hirsch, M. Brettreich, Fullerenes: Chemistry and Reactions, Wiley-VCH, (2005)
3. J. Oschwald, M. B. Minameyer, K. L. Maxouti, A. Hirsch, T. Drewello, submitted.

WP-172 / GAS PHASE STRUCTURES AND THERMOCHEMICAL PARAMETERS OF PROTONATED 5-HMF ISOMERS

Federico Pepi⁽¹⁾ - *Stefania Garzoli*⁽¹⁾ - *Pierluigi Giacomello*⁽¹⁾ - *Giulia de Petris*⁽¹⁾ - *Anna Troiani*⁽¹⁾ - *Rino Ragno*⁽¹⁾ - *Alexandros Patsilinos*⁽¹⁾ - *Lorenzo Antonini*⁽¹⁾ - *Andreina Ricci*⁽²⁾ - *Chiara Salvitti*⁽¹⁾

"Sapienza" University of Rome, Department of Chemistry and Drugs Technologies, Rome (1) - University of Campania, yet Second University of Naples, Department of Math and Phys, Caserta (2)

Keywords: Platform molecules, Biomass, 5-HMF, proton affinity, gas phase ion chemistry

Introduction

5-hydroxymethyl-2-furaldehyde (5-HMF) is produced in high yields from the acid catalyzed dehydration of carbohydrates. This reaction represents the key process of the cellulosic biomass conversion into furan-type platform chemicals. [1] Recently, we have employed mass spectrometric techniques to investigate the acid-catalyzed D-glucose and D-fructose dehydration highlighting the formation in the gas-phase of a protonated 5-HMF isomeric mixed population[2-3].

Methods

Full scan and MS_n mass spectra were acquired using a Bruker Amazon SL ion trap (IT) and a Waters Ultima

Q-TOF operating in the positive ion mode. Proton affinity (PA) and gas-phase basicity (GB) of 5-HMF were obtained by using the Wesdemiotis and Fenselau Cooks's "extended" kinetic method. The PA and GB values as well as the geometries of protonated isomeric structures were fully optimized in the gas-phase at the B3LYP/6-311++G** level of theory.

Results

In order to clarify the stability and structures of protonated 5-HMF isomers, the geometries of six different protomers were optimized at the B3LYP/6-311++G** level of theory and the energetics of their dissociation compared with the experimentally observed mass spectrometric fragmentation. Moreover, the unknown PA and GB of 5-HMF were measured by the extended Cooks's kinetic method along with theoretical calculations. The experimental PA value of 207.2 ± 3 kcal mol⁻¹ is in excellent agreement with the value of 207.5 kcal mol⁻¹ computed at the B3LYP/6-311++G** level of theory. The whole picture emerging from experimental and theoretical results allowed to attribute the structures of different protonated 5-HMF isomers to the *m/z* 127 ions generated by the acid catalyzed decomposition of hexose carbohydrates in the gas phase.

Conclusions

The unknown proton affinity and gas phase basicity of 5-HMF were determined by the joint application of Cook's kinetic method and quantum mechanical calculations. Theoretical calculations identify the oxygen atom of the aldehydic group as the most basic site. The formation of less stable protonated 5-HMF isomers as final product of the acid catalyzed decomposition of D-fructose and D-glucose was demonstrated on the basis of energetic considerations.

Novel Aspect

In this work the unexplored gas-phase ion chemistry of protonated 5-HMF, one of the top ten bio-based platform chemicals, was firstly investigated.

References

1. Ragauskas A. J., Williams C. K., Davison B. H., Britovsek G., Cairney J., Eckert C. A., Frederick W. J. Jr., Hallett J. P., Leak D. J., Liotta C. L., Mielenz J. R., Murphy R., Templer R., Tschaplinski T., *Science* 311, 484-489 (2006).
2. Ricci A., Di Rienzo B., Pepi F., Troiani A., Garzoli S., Giacomello P., *J. Mass Spectrom.* 50, 228–234 (2015).
3. Pepi F., Ricci A., Garzoli S., Troiani A., Salvitti C., Di Rienzo B., Giacomello P., *Carbohydr. Res.* 413, 145-150 (2015).

WP-173 / SELECTIVE CONVERSION OF PROTONATED D-FRUCTOSE TO 5-HMF AND 2-FA: A NITROGEN-BASE ASSISTED DEHYDRATION PROCESS IN THE GAS-PHASE

Anna Troiani⁽¹⁾ - **Giulia de Petris**⁽¹⁾ - **Federico Pepi**⁽¹⁾ - **Stefania Garzoli**⁽¹⁾ - **Marzio Rosi**⁽²⁾ - **Andreina Ricci**⁽³⁾ - **Chiara Salvitti**⁽¹⁾

Università di Roma " Sapienza", Dip. di Chimica e Tecnologie del Farmaco, Rome (1) - Dip di Ingegneria Ambientale e Civile, Università di Perugia, Perugia (2) - Università della Campania "L. Vanvitelli", Dip. di Matematica e Fisica, Caserta (3)

Keywords: fructose dehydration, Biomass, 5-HMF, gas phase ion chemistry, computational **Methods**

Introduction

Understanding on a molecular level, the dehydration mechanism of sugar monomers from hemicellulose and cellulose, (e. g xylose, glucose, fructose) the main constituents of lignocellulosic biomass, to platform molecules such as 5-hydroxymethyl 2-furaldehyde (5- HMF) and 2-furaldehyde (2- FA), can offer the possibility of increasing selectivity and reaction yields in solution for the exploitation of this sustainable source.

Methods

A joined approach based on mass spectrometric techniques (Electrospray Triple Quadrupole Mass Spectrometry (ESI-TQ/MS) and Electrospray Ion Trap Mass Spectrometry (ESI-QIT/MS)) together with theoretical calculations at the B3LYP/6-31 + G(d,p) level of theory was used for the study of the gas phase selective conversion of protonated D-fructose to 5-HMF and to 2-FA

Results

Our mass spectrometric approach for the gas phase study of acid catalyzed dehydration of xylose, fructose and glucose[1-3] allowed us to structurally characterize reaction intermediates and final products and unequivocally ascertain reaction sequences. Molecular level evidences provided by our mass spectrometric approach highlighted a gas phase pathway for the selective conversion of fructose to 5-HMF through a nitrogen assisted dehydration process.

Conclusions

Based on mass spectrometric evidences, theoretical calculations allowed a gas-phase mechanism for the selective conversion of protonated D-fructose to 5-HMF and 2-FA assisted by nitrogen-base, to be proposed.

Novel Aspect

In this work the gas phase selective fructose dehydration to 5-HMF and 2 FA assisted by nitrogen bases was firstly investigated.

References

1. Ricci A., Piccolella S., Pepi F., Garzoli S., Giacomello P., The mechanism of 2-furaldehyde formation from d-xylose dehydration in the gas phase . A tandem mass spectrometric study. J. Am. Soc. Mass Spectrom. 24, 1082-1089 (2013)
2. Ricci A., Di Rienzo B., Pepi F., Troiani A., Garzoli S., Giacomello P., Acid -catalyzed glucose dehydration in the gas phase: a mass spectrometric approach. J. Mass Spectrom. 50, 228–234 (2015)
3. Pepi F., Ricci A., Garzoli S., Troiani A., Salvitti C., Di Rienzo B., Giacomello P., A mass spectrometric study of the acid catalyzed d- fructose dehydration in the gas phase. Carbohydr. Res. 413, 145-150 (2015)

WP-174 / SULFUR DIOXIDE UPTAKE BY SODIUM CARBONATE CLUSTER ANIONS IN THE GAS-PHASE

Chiara Salvitti - Anna Troiani - Giulia de Petris

"Sapienza" University of Rome, Departement of Chemistry and Technology of Drugs, Rome

Keywords: ion-molecule reactions, gas-phase chemistry, sulfur dioxide, carbon dioxide, intra-cluster reaction

Introduction

Carbonate (CO_3^{2-}) and sulfate (SO_4^{2-}) species are reactive inorganic components of the fine particulate matter (PM) implicated in several environmental and human health issues. Sulfate aerosols usually arise from the atmospheric oxidation of SO_2 , contributing to acid rain and climate change. [1] In this study the reactivity of $[(\text{Na}_2\text{CO}_3)_n\text{NaCO}_3]^-$ cluster ions ($n \geq 1$) was probed towards SO_2 and $^{13}\text{CO}_2$, as a model of processes occurring at the liquid/gas phase interface.

Methods

The experiments were performed on a LTQ XL linear ion trap (Thermo Fisher Scientific) equipped with an electrospray ionization (ESI) source in negative ion mode. The instrument was in-house modified [2] to allow the introduction of neutral reagent gases (SO_2 , $^{13}\text{CO}_2$) into the ion trap and measure the kinetic rate

constants of the observed ion-molecule reactions. $[(\text{Na}_2\text{CO}_3)_n\text{NaCO}_3]^-$ cluster ions were generated by spraying a 10^{-3} M solution of Na_2CO_3 in $\text{H}_2\text{O}/\text{CH}_3\text{CN}$ 1:3.

Results

Cluster ions with the general formula $[(\text{Na}_2\text{CO}_3)_n\text{NaCO}_3]^-$ were observed in the gas-phase with $n \geq 1$. In the presence of SO_2 these carbonate species were quantitatively converted into the corresponding sulfite cluster ions $[(\text{Na}_2\text{SO}_3)_n\text{NaSO}_3]^-$ through a sequential replacement of each CO_2 moiety with a SO_2 molecule. The $\text{SO}_2 \rightarrow \text{CO}_2$ conversion may occur through the direct transfer of O_2^- from CO_3^{2-} to a SO_2 molecule, as reported in a previous study identifying the reactions of SO_2 at the surface of a molten carbonate eutectic. [3] The rate constants related to the SO_2 adsorption and CO_2 release were measured by monitoring the signal of the selected carbonate cluster ion as a function of SO_2 concentration. The obtained values were investigated on the basis of n . All the reactions are very fast and efficient. The same experiments were performed also in the presence of $^{13}\text{CO}_2$. As in the case of SO_2 , $^{13}\text{CO}_2$ is incorporated into the cluster ion structures by replacing CO_2 and leading to $[(\text{Na}_2^{13}\text{CO}_3)_n\text{Na}^{13}\text{CO}_3]^-$ ions with rate constants lower than those obtained for SO_2 .

Conclusions

SO_2 molecules were converted into gaseous CO_2 via the reaction of carbonate cluster ions and the consequent formation of sulfite cluster species. Although varying the number n of (Na_2CO_3) moiety of the reactant ion, CO_2 is always efficiently replaced and SO_2 entrapped. The same reactions were observed also in the presence of labeled $^{13}\text{CO}_2$, but showing rate constants lower than those measured for SO_2 .

Novel Aspect

A liquid/gas phase model has been reported to describe the processes involved in the maintenance of the atmospheric SO_2/CO_2 balance or in the SO_2 removal from flue gases.

References

1. Zhang R., et al., Chem. Rev. 115, 3803 (2015).
2. de Petris G., Cartoni A., Troiani A., Barone V., Cimino P., Angelini G., Ursini O., Chem. Eur. J. 16, 6234 (2010).
3. Krebs T., Nathanson G. M., Proc. Natl. Acad. Sci. 107, 6622 (2010).

WP-175 / DEVELOPMENT OF TUNER INTEGRATED MICROWAVE DRIVEN RADICAL BEAM SOURCES FOR RADICAL INDUCED DISSOCIATION MOLECULAR STRUCTURAL ANALYSIS

Yuji Shimabukuro⁽¹⁾ - ***Hidenori Takahashi***⁽²⁾ - ***Shinichi Iwamoto***⁽²⁾ - ***Koichi Tanaka***⁽²⁾ - ***Motoi Wada***⁽¹⁾

Plasma Physics Laboratory, Doshisha University, Kyotanabe (1) - Koichi Tanaka Mass Spectrometry Research Laboratory, Shimadzu Corporation, Kyoto (2)

Keywords: Microwave Plasma, Radical Induced Dissociation, Tandem Mass Spectrometer, Post Translational Modification, Phospholipid

Introduction

We have been developing a new mass analysis system by coupling a tandem mass spectrometer to a radical beam source. This research is an extension of our previous work: Hydrogen Attachment/Abstraction Dissociation (HAD) for peptide fragmentation[1]. A wide variety of radicals produced in the microwave plasma is introduced to a QIT confining peptide and phospholipid ions to investigate the possibilities of the gas phase reactions for biomolecular analyses.

Methods

Two types of tuner integrated microwave driven radical sources (CCP [2] and ECR-ICP) have been developed

to produce electrically neutral species from any reactive gases such as oxygen or water vapor. These radical sources directly deliver radical beams to a QIT chamber of MALDI QIT-TOF mass spectrometer to investigate the reaction between trapped analyte ions and injected species. In this report, substance P and phospholipids were used as analyte ions.

Results

Both radical sources can commence hydrogen, nitrogen, oxygen, methane, and water discharges. The CCP source forms the plasma plume which realizes high transport efficiency under a higher gas pressure and high input power condition. Meanwhile the ICP source can ignite and sustain the discharge at a lower pressure and lower input power. The OES results of hydrogen ICP exhibits the high degree of dissociation: prominent atomic lines with faint molecular band spectra. The HAD reaction rates for substance P with hydrogen radicals generated by three types of radical sources are investigated to evaluate the radical source performance. Hydrogen radicals formed in the pure hydrogen and water vapor discharge of CCP source can realize the HAD cleavage similar to the one obtained by a thermal cracking source, while the ICP source did not cleave the precursor ions. On the other hand, when the target ions are phospholipids, the injection of reactive species from the water vapor discharge from both CCP and ICP can realize the double-bond specific fragmentations.

Conclusions

H attachment pattern for QIT confined C_{60} ions are found different among thermal cracking, CCP and ICP sources. A series of experiments has shown that the CCP source produces a plasma of a higher effective temperature appropriate for peptide analysis that requires higher reaction energy. Meanwhile the ICP source that generates lower effective temperature plasma enables phospholipids analysis with the low energy radicals for double bond specific fragmentation.

Novel Aspect

The versatile radical sources which have wide range of energy will open new approaches of radical attachment/abduction mass spectrometry.

References

1. Takahashi H., Sekiya S., Nishikaze T., Kodera K., Iwamoto S., Wada M., Tanaka K., Anal. Chem. 88, 3810 (2016).

WP-176 / CONTROLLED REACTIONS OF METHYL CATION WITH OXYGEN COMPOUNDS AND HYDROCARBONS

Roland Thissen - Claire Romanzin - Christian Alcaraz

Laboratoire Chimie Physique, Université Paris Saclay, Orsay

Keywords: Ion-molecule reaction, Photoionization, Planetary science, Absolute cross section, Physical Chemistry

Introduction

The methyl carbocation is a reactive species whose presence has been detected in gaseous environments fed by high energy sources, such as the interstellar medium, planetary ionospheres and laboratory plasmas. Its reactivity is described only in its ground state while excited species are present and might affect the chemistry. Our established method of photoionization allows to study and quantify effects of internal and collision energy on reactivity.

Methods

The CERISES setup (a Quad-Oct-Quad guided ion-beam experiment) allows studying the effects of internal

energy (vibrational, electronic) on the reactivity of cations with neutral gaseous species. Ions with controlled degree of internal excitation are produced by photoionization at the DESIRS beamline of the SOLEIL synchrotron [2,3]. After mass selection, reaction occurs in calibrated pressure of target and products are then analysed in mass and counted.

Results

The $\text{CH}_3\cdot$ radicals are formed by pyrolysis in a molecular beam of the nitromethane (CH_3NO_2) precursor molecule [1]. The CH_3^+ ion is produced via direct photoionization of the $\text{CH}_3\cdot$ radicals. This allows to define the average internal energy of the ion. Hence, we have studied the spectroscopy of the cation in order to derive the effective internal energy of the cation at all available photon energies (from 9.8 to 15.5 eV). In a second step, the reactivity of CH_3^+ with (i) formic acid HC(O)OH (ii) methanol CH_3OH , and (iii) 10 hydrocarbon molecules (ranging from C1 to C4 with isomers and different levels of unsaturation) has been studied. When necessary, isotopomers were used to clarify reaction channels. Branching ratios and absolute reaction cross-sections for each systems and channel are obtained and the dependence on the vibrational and electronic excitation as well as kinetic energy (0-20 eV, lab) of the parent ion are described. Absolute reaction cross sections are obtained for all the systems allowing to derive a coherent image of the reaction dynamics.

Conclusions

Strong effects of internal energy (vibrational as well as electronic) are observed in all the considered systems. New reaction channels appear and branching ratios are perturbed, allowing an improved description of the actual chemistry that occurs in terrestrial or astrophysical plasma. Results are analyzed for each system in terms of partition between charge transfer, dissociative charge transfer, and new bond formation.

Novel Aspect

Access to direct photoionization of radicals, species which are important and difficult to handle, to probe effects of internal energy in ion molecule reactivity.

References

1. B. Cunha de Miranda, C. Romanzin, S. Chefdeville, V. Vuitton, J. Žabka, M. Polášek, C. Alcaraz, J. Phys. Chem. A 119, 6082 (2015).
2. L. Nahon et al., J. Synchrotron Radiat. 19, 508 (2012).
3. A. Cernuto, A. Lopes, C. Romanzin, B. Cunha de Miranda, D. Ascenzi, P. Tosi, G. Tonachini, A. Maranzana, M. Polasek, J. Zabka, C. Alcaraz, J. Chem. Phys. 147, 154302 (2017).

WP-177 / TWO MOLECULES ARE BETTER THAN ONE: COOPERATIVE EFFECTS OF TWO SULFUR DIOXIDE AND SULFUR DIOXIDE/WATER IN BOND FORMING AND HYDROLYSIS REACTIONS OF POLYVANADATE DIANIONS IN THE GAS PHASE

Anna Troiani⁽¹⁾ - **Marzio Rosi**⁽²⁾ - **Stefania Garzoli**⁽¹⁾ - **Chiara Salvitti**⁽¹⁾ - **Giulia de Petris**⁽¹⁾

“Sapienza” University of Rome, Dipartimento di Chimica e Tecnologie del Farmaco, Rome (1) - University of Perugia, Dipartimento di Ingegneria Civile e Ambientale and ISTM-CNR, Perugia (2)

Keywords: Polyvanadate doubly-charged anions, sulfur dioxide, bond-forming reactions, hydrolysis, cooperative effects

Introduction

Vanadium oxides are mostly used as catalysts in many oxidation reactions such as of SO_2 to SO_3 , whereas hydrated vanadium oxides have potential application in Li-ion batteries [1, 2]. Following our interest in the

chemistry of SO₂, we recently reported an unprecedented bond-forming and hydrolysis reactions of doubly-charged vanadium oxide and hydroxide cluster anions promoted by SO₂, leading to singly-charged ions, as a result of a charge separation process [3, 4].

Methods

The reactions have been studied in the gas phase at room temperature by an in-house modified linear quadrupole ion trap mass spectrometer. NaV₃O₉²⁻ and H₂NaV₃O₁₀²⁻ dianions, prepared by electrospray ionization (ESI) of NaVO₃ aqueous solutions, were reacted in the gas-phase with SO₂. The rate constant, branching ratio and the efficiency of the reactions have been measured.

Results

NaV₃O₉²⁻ reacts with sulfur dioxide giving a first stable association product [NaV₃O₉SO₂]²⁻, that further reacts with a second molecule of SO₂ forming two singly charged products, NaV₂O₆SO₂⁻ and VO₃SO₂⁻. The reaction involves formation of new V-O and S-O bonds, by breaking a V-O bond containing a bridging oxygen atom. By contrast, H₂NaV₃O₁₀²⁻, the hydrated form of NaV₃O₉²⁻, reacts with only one SO₂ also giving the hydrolysis products HOSO₂⁻ and HNaV₃O₉⁻ in addition to those arising from the bond-forming process, H₂VO₄⁻ and NaV₂O₆SO₂⁻. Both reactions proceed through the formation of long-lived intermediates, where the bonding to the neutral SO₂ is favoured with respect to the charge separation, that would arise from an electron transfer to SO₂. Instead, the charge separation is the consequence of a chemical reaction leading to the formation of two singly charged products. A common feature is the crucial role played by sulfur dioxide in promoting charge separation processes.

Conclusions

The sequential addition of two molecules of sulfur dioxide to thermal doubly-charged vanadium oxide anions, NaV₃O₉²⁻, results in the breaking of the stable V₃O₉ bone, whereas in the case of the hydroxide H₂NaV₃O₁₀²⁻ anions one SO₂ molecule is sufficient to promote also hydrolysis. These fast and efficient reactions are associated with charge separation processes that result in two singly charged ions by formation of new V-O and S-O bonds.

Novel Aspect

Cooperativity of two SO₂ or SO₂/H₂O molecules in the bonding to polyvanadate dianions has proved fundamental in addressing bond-forming or hydrolysis reactions.

References

1. Weckhuysen B. M., Keller D. E., Catal. Today78, 25 (2003)
2. Zhang C., Song H., Zhang C., Liu C., Liu Y., Cao G., J. Phys. Chem. C 119, 11391 (2015)
3. Troiani A., Rosi M. Garzoli S., Salvitti C., de Petris G., Chem. Eur. J. 23, 11752 (2017)
4. Troiani A., Rosi M. Garzoli S., Salvitti C., de Petris G. New J. Chem., 42, 4008 (2018)

WP-178 / GAS-PHASE FRAGMENTATION MECHANISMS STUDIES OF DOUBLY-PROTONATED PEPTIDE IONS CONTAINING BASIC AMINO ACID RESIDUES BY MASS SPECTROMETRY

Ahmet Emin Atik⁽¹⁾ - **Talat Yalcin**⁽²⁾

Turgut Pharmaceuticals, Turgut Pharmaceuticals, Istanbul (1) - IzTech., Chemistry, Urla-Izmir (2)

Keywords: Basic Amino Acid Residues, Doubly-Protonated Peptide Ions, Gas Phase, Peptide Fragmentation

Introduction

In recent years, many studies have been carried out for the gas-phase reaction mechanisms of the protonated peptides using mass spectrometry and many of these studies are related with singly-

protonated peptides (1). However, there are a few studies related to multiply-protonated peptides in the literature (2,3). Most of these studies are not systematic studies, therefore they are not sufficient enough to better understand the gas-phase reaction mechanisms of peptides. Multiply-protonated peptides can be generated in ESI-MS studies if the peptides contain a basic amino acid residues.

In order to determine the amino acid sequence of the peptides correctly and reliable, the gas-phase reaction mechanism of peptides should be studied and understood very well. Multiply-protonated peptides can be easily formed if the peptide sequence consists of basic amino acid residue. As all know, when the proteins undergo enzymatic cleavage (especially with Trypsin), the C-terminal of truncated peptides always contain a basic amino acids such as lysine (K) or arginine (R). Especially, the gas-phase fragmentation mechanisms of multiply-protonated peptides should be investigated in details by MS to elucidate the sequence as well as the function of proteins.

Methods

Mass spectrometry experiments were conducted on a LTQ XL linear ion-trap mass spectrometer (Thermo Finnigan, San Jose, CA, USA) equipped with an electrospray ionization (ESI) source. A 100 pmol L⁻¹ peptide solutions were introduced into the ion source at a flow rate of 5 L/min. The peptide working solutions were prepared in MeOH/dH₂O/HCOOH (50:50:1, v/v/v) containing mixture. Full scan MS spectra were acquired in the *m/z* range 150–1000.

Results

During the study, XXXXXXH-NH₂ (X = A, C, F, G, L, M, S, T, V, Y) and XYAGFLV-NH₂, YXAGFLV-NH₂, YAXGFLV-NH₂, YAGXFLV-NH₂, YAGFXLV-NH₂, YAGFLVX-NH₂ peptide derivatives were studied. In literature, it has been reported that peptide chain length can effect the fragmentation mechanisms of doubly-protonated peptides using AlaxH (X= 5,6,7,8 and 10).

In this work, the peptide chain length is kept constant and the position of basic residues, lysine (K), arginine (R), and histidine (H), is changed within the peptide sequence to probe the effects of position as well as the identity of these amino acid residues on the fragmentation pathways. The effects of location of basic residues (close to the C-terminal and N-terminal) was also studied.

Conclusions

This work presents a detailed investigation of multiply charge ions produced different peptide series. It was shown that, all the peptide series were showed the class 1 ions (*b*₂/*y*₅) considering basic amino acid position starting from N-position (3,4,5,6,and 7) *b*₂ and *y*₅ ion intensities were changed with respect to position of basic amino acid. The formation of *b*₂ and *y*₅ ion was observed at 4th position at the peptide series of XYAGFLV-NH₂, YXAGFLV-NH₂, YAXGFLV-NH₂, YAGXFLV-NH₂, YAGFXLV-NH₂, and YAGFLVX-NH₂ (where X= K, H and R) not the 3rd position like the peptide series of XAAAAAA-NH₂, AXAAAAA-NH₂, AAXAAAA-NH₂, AAAAXAA-NH₂, AAAAXA-NH₂, and AAAAAAX-NH₂ (where X= K, H, and R) *b*₂/*y*₅ (class 1) ion formation was affected by peptide sequences. Especially, the ions (*b*₂/*y*₅) produced from SSSSSH-NH₂.

Novel Aspect

The purpose of the current study is the investigation of the gas-phase fragmentation mechanisms of multiply-protonated peptides to gain peptide sequence correct and reliable.

References

1. Tang, X. J., Thibault, P., Boyd, R. K. 1993, Analytical Chemistry, 65, 2824-2834, (1993)
2. Savitski, M. M., Fälth, M., Fung, Y. M. E., Adams, C. M., Zubarev, R. A., Journal of the American Society for Mass Spectrometry, 19, 1755-1763, (2008)
3. Harrison, A. G., Journal of the American Society for Mass Spectrometry, 20, 1890-1895, (2009)

Acknowledgements

The authors gratefully acknowledge to Scientific and Technological Research Council of Turkey TUBITAK, for financial support under project 215Z257 and also grateful the State Planning Organization (DPT) of Turkey

for the funding of the Biological Mass Spectrometry and Proteomics Facility located at Izmir Institute of Technology.

WPS-S08 - FTMS and innovative methods for (bio)analyses

Chairs: Ljiljana Paša-Tolić, Giuliana Bianco

WP-179 / TOP-DOWN ANALYSIS OF ENDOGENOUS TAU FRAGMENTS IN HUMAN CEREBROSPINAL FLUID USING IMMUNOPRECIPITATION AND FOURIER TRANSFORM MASS SPECTROMETRY

Gunnar Brinkmalm - Claudia Cicognola - Johan Gobom - Erik Portelius - Henrik Zetterberg - Kaj Blennow - K. Höglund

University of Gothenburg, Sahlgrenska Academy, Mölndal, Sweden

Keywords: top-down MS, immunoprecipitation, tau protein, cerebrospinal fluid, endogenous peptides

Introduction

Many proteins in cerebrospinal fluid (CSF) exist in processed forms, many of which reflect corresponding species present in the brain. Previously an abundance of endogenous fragments from amyloid precursor protein, including variants of amyloid β has been identified both in CSF and brain. Until recently, the other core biomarker for Alzheimer's disease, tau, has been more difficult to assess with mass spectrometry (MS) due to its lower abundance in CSF.

Methods

Tau was enriched from human CSF by immunoprecipitation (IP) using six different antibodies directed towards different tau regions. Top-down was the main approach, but bottom-up was also utilised for these samples. Analysis was performed by nanoflow liquid chromatography–Fourier transform MS. Quantitative data were obtained by addition of an isotope labeled protein standard before IP. Data processing was performed using Mascot Distiller and PEAKS Studio.

Results

Starting with three commercially available antibodies endogenous tau fragments were identified using top-down MS of IP'd CSF samples. Based on these data two new neospecific antibodies (Tau123_C and Tau224_C) were developed. In total more than 100 endogenous tau proteoforms were identified by IP-MS using all six antibodies. These mainly covered the N-terminal half, but there were also species originating from the C-terminal quarter. The protein N-terminal was found to be acetylated and phosphorylation at Thr181 was observed. The bottom-up approach tryptic peptides covering slightly more than the N-terminal half were identified. From the quantitative measurements it could be observed that proteoforms covering the mid-region are most abundant, followed by the N-terminal region and a small contribution from the C-terminal part. Three major cleavage sites were identified; in the N-terminal quarter, in the mid-region, and in the C-terminal quarter. These results are in concordance with published data indicating that CSF tau appears in various processed forms, dominated by species originating from the N-terminal half of the protein [1].

Conclusions

IP-MS is a powerful tool for proteoform mapping in biological samples. It is vital to understand tau processing and its relation to neurodegenerative diseases. These data can be used to guide development of

new immunoassays targeting specific tau proteoforms.

Novel Aspect

Characterisation of endogenous tau fragments in CSF using IP-MS methodology reveals the complexity of biomarker development and evaluation.

References

1. NR Barthélemy, F Fenaille, C Hirtz, N Sergeant, S Schraen-Maschke, J Vialaret, L Buée, A Gabelle, C Junot, S Lehmann, F Becher, *J Proteome Res* 15, 667-676(2016).

WP-180 / EXTREME RESOLUTION FT-ICR-MS TOP-DOWN ANALYSIS OF GLYCATED PROTEINS

Carlos Cordeiro

Laboratório de FT-ICR e Espectrometria de Massa Estrutural, Faculdade de Ciências, Universidade de Lisboa, Lisboa

Keywords: FT-ICR-MS, Protein Glycation, ECD, Absorption mode

Introduction

Methylglyoxal, a non-enzymatic by-product of glycolysis, is the most relevant glycation agent in vivo. Protein modification by methylglyoxal is found in association to several human pathologies, most notably amyloid neurodegenerative diseases [1]. Characterization of glycated proteins is extremely challenging [2] calling for extreme resolving power and multiple ion-fragmentation methods.

Methods

Methylglyoxal modified proteins were analysed by direct infusion in a 7 Tesla Bruker Solarix XR mass spectrometer equipped with the dynamically harmonized FT-ICR ParaCell, in magnitude and absorption modes, at 4 to 8 MB acquisition sizes. Proteins were fragmented by using different ion-activation methods, including ECD.

Results

Top down analysis of glycated proteins allows the identification of most modified species at MS, while MSMS analysis pinpoints location and chemical nature of modifications. While in vitro glycated proteins are highly heterogeneous, in vivo glycated ones show a higher degree of homogeneity. Most of the methylglyoxal derived glycation products can be identified and their location defined.

Conclusions

Glycated proteins have often been analysed by bottom-up methods that do not allow any insight on the diversity and heterogeneity. Top-down analysis reveals the extreme complexity of glycated proteoforms, paving the way for interpreting the glycation effects on protein structure and function, essential for understanding its effects in living cells.

Novel Aspect

Extreme resolution is a key requirement for analysing heterogeneously modified proteins, combined with a diversity of MSMS methods and in cell isolation of modified proteoforms.

References

1. Vicente Miranda, H et al., *Brain*, 1401399-1419 (2017)
2. Gomes, RA et al., *Biochem J*, 416 317-326 (2008).

WP-181 / IDENTIFICATION OF NEW LIGAND-PROTEIN COMPLEX FOR MYCOBACTERIUM TUBERCULOSIS USING FOURIER TRANSFORM MASS SPECTROMETRY

Angela Di Capua⁽¹⁾ - ***Peter Myler***⁽²⁾ - ***Ronald Quinn***⁽¹⁾

Griffith Institute for Drug Discovery, Don Young Road, Brisbane (1) - Center for Infectious Disease Research, University of Washington, Seattle (2)

Keywords: natural products, native mass spectrometry, target identification

Introduction

Natural products are an infinite resource for drug discovery and development. Native mass spectrometry (MS) is an approach used to study proteins in the native folded state and useful to investigate non-covalent ligand protein complexes. We applied this method to screen our natural-product-based fragment library against Malaria proteins [1] and natural product fraction library against Tuberculosis (TB) proteins.

Methods

37 frozen, shipped proteins are buffer-exchanged with ammonium acetate solution and then diluted to run the mass spectrometry assay. The concentration range used for protein is between 1-10 μ M. Screening compounds or fractions are pooled together, dried and treated with 1 μ L methanol. The molar ratio of protein to single compound is 1:5. Screening is performed using Triversa Nanomate (Advion) coupled with a FT-ICR-MS Solarix (Bruker).

Results

37 proteins were tested for visibility in native state using FT-ICR-MS. 31 proteins have been screened against 362 pooled active fractions. 26 proteins are able to bind a component with high (ligand-protein complex > 50%), medium (ligand-protein complex > 15%), or weak (ligand-protein complex < 15%) affinity. We identified components able to bind only one protein, components able to bind different proteins and proteins that show different ligand-protein complexes from the same pool fraction. Two proteins showed high ligand ability (more than 7 binding complexes), 16 proteins showed medium ligand ability (more than 3 binding complexes) and 7 proteins showed low ligand ability (less than 3 binding complex). Fraction were confirmed and compound isolation is in progress.

Conclusions

We used the native mass spectrometry approach for a high throughput screening with natural fractions active against TB. We investigate 37 potential targets and we have identified 26 proteins that bind to a total of 90 components from active tuberculosis fractions.

Novel Aspect

Native mass spectrometry can be considered an alternative approach for High Throughput Screening of pure compounds as well as fractions libraries and target validation.

References

1. Vu, H., et al., ACS Infect. Dis. 2018, 4, 431-444

WP-182 / COUPLING GAS-PHASE FRACTIONATION AND PROTON-TRANSFER REACTION FOR IMPROVED TOP-DOWN PROTEOMICS

Luca Fornelli⁽¹⁾ - ***Romain Huguet***⁽²⁾ - ***Christopher Mullen***⁽³⁾ - ***Ryan Fellers***⁽¹⁾ - ***Rafael Melani***⁽¹⁾ - ***John Syka***⁽³⁾ - ***Neil Kelleher***⁽¹⁾

Northwestern University, Molecular biosciences/Chemistry, Evanston (1) - Thermo Scientific, LSMS Product

Management, San Jose (2) - Thermo Scientific, LSMS R&D, San Jose (3)

Keywords: proton-transfer reaction; top-down; proteomics; Orbitrap; gas-phase fractionation.

Introduction

The major limitation in top-down proteomics is its reduced throughput. Here, we implement a series of gas-phase fractionation strategies on a Thermo Scientific Orbitrap Fusion Lumos Tribrid mass spectrometer to improve the number of detectable proteoforms. Specifically, we applied proton-transfer reaction (PTR) to both deepen the proteome coverage in the 0-30 kDa mass range and to facilitate the analysis of large proteoforms up to ~60 kDa.

Methods

Classic data-dependent acquisition (DDA) methods (*i.e.*, “top-N” strategies) based on either full MS1 or selected ion monitoring (SIM) scans for precursor identification were compared with gas-phase fractionation based on PTR. Quadrupole-isolated (width: 10-50 m/z units) multiply charged protein ions were subjected to PTR MS2 in the ion trap; precursors identified by PTR MS2 were isolated for HCD MS3 fragmentation using a MultiNotch method in the linear trap.

Results

Orbitrap MS allows for the characterization of proteoforms from 0 to ~60 kDa. [1] However, even when applying extensive fractionation only the most abundant proteoforms are typically characterized. We performed gas-phase fractionation as previously described by the McLuckey and Hunt groups, by applying PTR on a selected portion (10-50 m/z units) of the spectrum at a time, to reduce spectral complexity. By reducing the charge state of protein cations, PTR also limits the problems in the narrow isolation of single proteoform precursors ahead of dissociation. We analyzed both bacterial (*Pseudomonas aeruginosa*) and human (HeLa cells) proteomes, in order to determine the different effects of PTR-based gas-phase fractionation on samples with different complexity. We also modified the synchronous precursor selection (SPS) method to adapt it to intact protein cations of different masses and charge states, making the selection of multiple precursors combined with HCD fragmentation efficient even for large proteoforms (>30 kDa).

Conclusions

We implemented a novel data acquisition strategy based on PTR MS2 followed by isolation of multiple precursors for MS3 fragmentation that is aimed at improving the proteome coverage in top-down proteomics. The flexibility offered by the Tribrid Orbitrap platform offers the possibility of adjusting the degree of gas-phase fractionation in consideration of the complexity of the targeted proteome.

Novel Aspect

Proton-transfer reaction is performed on a Tribrid Orbitrap mass spectrometer to detect and characterize more proteoforms in top-down proteomics.

References

1. Fornelli L., Durbin K. R., Fellers R. T., Early B. P., Greer J. B., LeDuc R. D., Compton P. D., Kelleher N. L., Journal of Proteome Research 16, 609-618 (2017).

WP-183 / DETAILED ANALYSIS OF TEMPORAL ION CLOUD DISTRIBUTIONS UNDERGOING COULOMB INTERACTIONS IN FT-ICR MASS SPECTROMETRY: A MANY-PARTICLE SIMULATION

Makoto Fujiwara - Hisashi Takahashi - Naohisa Happo - Koichi Tanaka

Hiroshima City University, Grad School of Information Sciences, Hiroshima

Keywords: FT-ICR, Coulomb Interaction, Many-Particle Simulation

Introduction

In FT-ICR MS, many different m/z ions are simultaneously trapped and excited in the cell and the temporal variation of each m/z ion cloud distribution is complicated. In particular, at high ion density the Coulomb interactions are important. To analyze the Coulomb-interaction-induced effects by many-particle simulation, we have combined a special-purpose computer for many-body problems with a PC and the computer configuration has been accelerated [1].

Methods

On the above accelerated configuration, the influence of Coulomb-interaction-induced degradation of excited coherent ion cyclotron motion on FT-ICR mass spectral peak shape was investigated by many-particle simulation. The simulation particles were equally assigned to each of m/z 130.0 and 130.1 ions. Ion excitation was performed by frequency sweep at 7 T and the mass spectra were simulated as a function of total ion population initially trapped in the cell.

Results

From analysis for temporal variations of the ion cloud distributions during detection period it has been found that the coherent ion cyclotron motion of m/z 130.0 ion (lighter weight ion) cloud is out of phase earlier than that of m/z 130.1 ion cloud with increasing initial ion population. However, under certain conditions it has also been found that the spectral peak shape obtained from m/z 130.1 ion cloud appears to be distorted earlier than that from m/z 130.0 ion cloud. We have extracted respective ICR signal components resulting from m/z 130.0 and 130.1 ion clouds. In comparison between the two components, it has been found that the component from m/z 130.1 ion cloud is more significantly modulated than that from m/z 130.0 ion cloud. Further we have analyzed in more precise detail temporal variations of the ion cloud distributions and have presumed that the phase locking may occur for a part of remained (non-ejected) m/z 130.1 ion cloud with m/z 130.0 ion cloud.

Conclusions

The simulation results suggest that the relationship between the Coulomb-interaction-induced degradation and the spectral peak shape is quite complicated.

Novel Aspect

Temporal variations of the ion cloud distributions undergoing Coulomb interactions in the FT-ICR cell have been analyzed in more detail by many-particle simulation.

References

1. M. Fujiwara, N. Happo, K. Tanaka, J. Mass Spectrom. Soc. Jpn. 58, 169-173 (2010).

WP-184 / FT-ICR MS WITH ABSORPTION-MODE DATA ACQUISITION: METHOD DEVELOPMENT AND HIGH-RESOLUTION APPLICATIONS

Anton N. Kozhinov - Konstantin O. Nagornov - Yury O. Tsybin

Spectroswiss Sàrl, EPFL Innovation Park, Lausanne

Keywords: FT-ICR MS, transient, absorption mode FT, data acquisition, ion phase

Introduction

Unlike the modern Orbitraps, many FT-ICR instruments worldwide offer magnitude-mode FT as a default, and sometimes the only, method to generate mass spectra. Previously, we introduced a data acquisition

technology that enables recording FT-ICR mass spectra directly in absorption mode, and tested this with a 7T LTQFT Ultra. Here, we extend this technology to more FT-ICRs, and evaluate it in wide-range applications including petroleomics, imaging, and proteomics.

Methods

Peptides, calibration mixtures, proteins, oils, and tissues were analyzed on commercial FT-ICR MS (9.4T solariX XR with ParaCell or custom NADEL cell, and two 7T LTQFT Ultra). Measurements were conducted in the w+ regime on the 7T and 9.4T instruments (original ICR cells), and in the wc and w+ regimes on the 9.4T instrument (NADEL cell). Transients were acquired in parallel with the instrument manufacturer's data, using a custom data acquisition system (FTMS Booster).

Results

A wide use of absorption-mode FT in FT-ICR MS applications is limited by high numerical costs of phase correction procedures, substantial phase distortions from standard data acquisition systems, and other complications. To avoid these, we developed a high-performance data acquisition system for recording FT-ICR data as transients with in-phase ions signals (or as mass spectra generated in absorption mode), compatible with both types of FT-ICRs from this work. The system operation was defined based on the textbook principles, namely: (i) the convolution property of FT and (ii) sampling of both detect and excite signals. The system was designed to sample and in-line process the two signals in order to produce a single transient for each measurement cycle. Our preliminary results demonstrate that with this technology FT-ICR mass spectra can be generated in absorption mode with no or a little tuning of the phase function. As a result, improved analyses of petroleum and proteomics samples, and increased throughput in tissue imaging are achieved.

Conclusions

This work's high-performance electronics and methods for acquiring FT-ICR mass spectra in absorption mode significantly simplify the use of absorption-mode FT processing in real-life high-resolution FT-ICR MS applications, reduce phasing artifacts in mass spectra, increase numbers of identified and quantified peptides, lipids and compound classes; and make for bringing the absorption mode processing to FT-ICR MS as a standard method to generate mass spectra.

Novel Aspect

Method development, implementation in hardware, and evaluation in applications of absorption-mode data acquisition for FT-ICR MS, including both, Bruker and Thermo, makes.

WP-185 / FACTOR XA ANTICOAGULANT BIOACTIVE PROTEINS DERIVED FROM SCORPION VENOM

Meng Li⁽¹⁾ - ***Pui Yiu Yuko Lam***⁽¹⁾ - ***Peng Chen***⁽²⁾ - ***Rémy Gavard***⁽³⁾ - ***Cookson K. C. Chiu***⁽¹⁾ - ***Christopher A. Wootton***⁽¹⁾ - ***Qiong Wu***⁽²⁾ - ***Mark P. Barrow***⁽¹⁾ - ***Hongzheng Fu***⁽²⁾ - ***Peter B. O'Connor***⁽¹⁾

University of Warwick, Department of Chemistry, Coventry (1) - Peking University, School of Pharmaceutical Science, Beijing (2) - University of Warwick, MAS CDT, Coventry (3)

Keywords: Fourier Transform Ion Cyclotron Resonance (FTICR) MS; Scorpion Venom; Anticoagulant drug; Proteomics; Chemometrics

Introduction

Of the 1400 species of scorpions worldwide [1], the Chinese scorpion is used in Chinese medicines for cardiovascular problems [2], antimicrobial [3], and tumors [4]. The venom is a complex mixture of proteins containing bioactive components [5], some showing inhibition of enzyme factor Xa, which prevents blood coagulation [6]. These proteins are often modified and crosslinked with multiple disulphide bonds, with

only partial genome sequences, therefore requiring de novo sequencing.

Methods

Crude scorpion venom was fractionated using step elution on C18 cartridge, and in-vitro inhibition of factor Xa was tested. Nano-LC-FTICR MS/MS analysis of each fraction was acquired with collisionally activated dissociation (CAD) and electron capture dissociation (ECD) as fragmentation methods. Chemometrics was used to obtain possible pharmaceutical targets. Spectra were collected on a 12 Tesla Bruker Solarix FTICR MS, using both top-down and bottom-up proteomics.

Results

Using direct infusion FT-ICR MS, more than hundred unique molecular species (approximately 2-12 kDa) were observed in the crude venom, but with bottom-up proteomics, only a small number of proteins can be identified, even though many of these unidentified peptides have good fragmentation spectra. To accurately identify potential bioactive compounds, methodologies were developed to correlate the spectra with the bioactivity results. Quantification of compounds present in each fraction made it possible to use linear regression approaches such as partial least square (PLS) to establish a correlation, which allowed for the identification of the compounds that have stronger correlations with the bioactivity. Therefore, a list of compounds that strongly correlate to the bioactivity could be potential pharmaceutical targets. In addition, from the current results, most of these target compounds are proteins that range from 2k to 5k Da, with 30% of the proteins have sequences that are not in the genome database from Uniport.

Conclusions

The potential factor Xa inhibitors were de novo sequenced using both top-down and bottom-up approaches. The high mass accuracy and resolving power of the FTICR-MS allowed for the highest possible confidence during analysis, which aids the sequencing study of these novel scorpion venom proteins. The method is applicable to identifying potential pharmaceutical leads from other natural products.

Novel Aspect

Ultra-high resolution FTICR MS for the analysis of scorpion venom reveals bioactive anticoagulant proteins.

References

1. Sridhara, S. ; Chakravarthy, A. K. ; Kalarani, V. ; Reddy, D. C., Diversity and Ecology of Scorpions: Evolutionary Success Through Venom, Springer Singapore, pp 57-80 (2016)
2. Wang, R. ; Moreau, P. ; Deschamps, A. ; de Champlain, J. ; Sauv e, R. ; Foucart, S. ; Bai, L. ; Lu, X. -R., Toxicon 32, 191-200 (1994).
3. Gao, B. ; Zhu S., Front. Microbiol. 9, 320 (2018).
4. D'Suze, G. ; Rosales, A. ; Salazar, V. ; Sevcik, C., Toxicon 56, 1497-1505 (2010).
5. Xu, X. ; Duan, Z. ; Di, Z. ; He, Y. ; Li, J. ; Li, Z. ; Xie, C. ; Zeng, X. ; Cao, Z. ; Wu, Y. ; Liang, S. ; Li, W., J. Proteom 106, 162-180 (2014).
6. Ortiz, E. ; Gurrola, G. B. ; Schwartz, E. F. ; Possani, L. D. Toxicon, 93, 125 (2015).
7. Xue, T. ; Ding, S. ; Guo, B. ; Zhou, Y. ; Sun, P. ; Wang, H. ; Chu, W. ; Gong, G. ; Wang, Y. ; Chen, X. ; Yang, Y., J. Med. Chem. 57, 7770-7791 (2014).

WP-186 / MICROCHIP CAPILLARY ELECTROPHORESIS-ESI-MS FOR RAPID, MULTI-LEVEL ANALYSIS OF COMPLEX PROTEINS

J. Scott Mellors - Erin Redman

908 Devices Inc., Life Sciences, Boston

Keywords: Protein characterization, capillary electrophoresis, microfluidics

Introduction

Mass spec characterization of proteins often requires multiple “levels” of analysis for complete characterization. For LC-MS methods these analyses typically require completely different columns and mobile phases. With Microchip CE-ESI-MS a single set of experimental conditions can be used for all levels. This enables rapid multilevel characterization of proteins, with no down-time between samples.

Methods

All work was performed using a commercially available microfluidic CE-ESI system (ZipChip, 908 Devices Inc.), attached to an orbitrap mass spectrometer (Thermo). Proteins were analyzed intact by diluting to an appropriate concentration before analysis. Middle-down and bottom-up analyses utilized standard reduction and enzymatic digestion methods. All samples were analyzed using the same microfluidic device and the same background electrolyte.

Results

To achieve successful analysis at all three structural levels, a BGE with a relatively low pH was chosen. In this denaturing BGE we observe just a single mAb peak, but the electrophoretic separation still functions to electrokinetically “desalt” the sample. This yields a clean mass spectrum with a simple dilute and shoot method. Middle-down analysis of mAbs yields separation of three main segments of the molecule: the free light chain (LC), and the Fc and Fd domains of the heavy chain. We also observe separation of variant forms of each of these segments. Top-down analysis of smaller proteins using this method yields similar separation of variants. The separation of sialic acid variants of alpha-1-acid glycoprotein will be presented as an example. Additionally, a high degree of sequence coverage in <20 minutes will be demonstrated for bottom-up analysis. results will focus on the ability to rapidly characterize protein structure and determine biologically relevant information.

Conclusions

The microfluidic CE-MS system is a simple and rapid way to perform multiple levels of protein characterization. By choosing a BGE compatible with all types of analysis, there is no need to change the system or analysis conditions between types of characterization assays. Therefore, users can obtain intact, middle-up, and peptide mapping data without any downtime due to system turnover.

Novel Aspect

Demonstration of a unique mass spec front end separation capability for rapid characterization of proteins

WP-187 / FLAVONOID ANNOTATION USING A PRODUCT ION-DEPENDENT MSⁿ DATA ACQUISITION METHOD ON A TRIBRID ORBITRAP MASS SPECTROMETER

Reiko Kiyonami⁽¹⁾ - *Iwao Sakane*⁽²⁾ - *Seema Sharma*⁽¹⁾ - *Graeme McAlister*⁽¹⁾ - *Amanda Souza*⁽¹⁾ - *Caroline Ding*⁽¹⁾ - *Andreas Huhmer*⁽¹⁾

Thermo Fisher Scientific, Chromatography and Mass Spectrometry, San Jose (1) - ITO EN, LTD, Central Research Institute, Tokyo (2)

Keywords: Flavonoid annotation; MSⁿ; Tribrid Orbitrap; Mass Frontier 8.0; Compounder Discoverer 3.0

Introduction

Flavonoid annotation from natural products remains challenging because of their structural diversity and lack of authentic standards. Simple MS² based analyses are often not sufficient for complete structural annotation of flavonoids. We present a new flavonoid annotation workflow that uses comprehensive fragment ion information from HCD MS-MS and higher order FTMSⁿ for rapid flavonoid annotation.

Methods

Flavonoid extracts from different types of natural products were analyzed using LC-MS. A modified Orbitrap Fusion™ Tribrid™ mass spectrometer was used for collecting HRAM MS and MSn data. The data were processed using Compound Discoverer™ 3.0 software. A novel structure ranking algorithm was applied to the MS and MSn data for confident structure elucidation of the unknown flavonoids.

Results

The flavonoid extracts from multiple tea and fruit/vegetable juice samples were analyzed using the developed flavonoid annotation workflow. The MSn (up to 5) scans were only triggered when one of the sugar neutral loss product ions (162.0523, 180.0628, 146.0574, 164.0685, 176.0315, and 194.0421) was detected in the MS/MS scan. The collected data were processed using Mass Frontier 8.0 and Compound Discoverer 3.0 software. A novel structure ranking algorithm included in the Compound Discoverer 3.0 software was applied to the MS-MS and MSn data for confident structure characterization of the unknown flavonoids based on ChemSpider database and custom flavonoids database. The MSn data were critical, especially for the annotation of flavonoid glycoconjugates and flavonoid isomers.

Conclusions

The new LC-MSn workflow allows intelligent MS/MS fragment ion dependent MSn data collection. In combination with dedicated software for data processing, the new workflow enables improved throughput, identification coverage and confidence for flavonoids annotation from natural products by using comprehensive fragment ion information from HCD MS-MS and higher order FTMSn.

Novel Aspect

The innovative LC-MSn workflow on a modified Tribrid instrument enables improved throughput, identification coverage and confidence for flavonoids profiling experiments.

WP-188 / FUNDAMENTALS OF FT-ICR MS AT THE CYCLOTRON FREQUENCY

Oleg Y. Tsybin⁽¹⁾ - ***Konstantin O. Nagornov***⁽²⁾ - ***Anton N. Kozhinov***⁽²⁾ - ***Yury O. Tsybin***⁽²⁾

Peter the Great Saint-Petersburg State Polytechnic University, Ion Physics Laboratory, Saint-Petersburg (1) - Spectroswiss, Biomolecular Mass Spectrometry Laboratory, Lausanne (2)

Keywords: FTMS, ion cyclotron resonance, ion motion, cyclotron frequency, transient

Introduction

Recently, we introduced and initially described FT-ICR MS at the cyclotron frequency instead of the reduced cyclotron frequency using the narrow aperture detection (NADEL) ICR cells and applicable for biomolecular applications [1-3]. Here, we will present an analytical model for signal generation at the cyclotron frequency based on our understanding of the underlying principles of this new concept of ion motion, excitation and detection in the NADEL ICR cells.

Methods

Transient waveforms at single and multiple cyclotron periods were measured experimentally and calculated theoretically by taking into account: (i) the trajectories of ion-oscillators with time-dependent instantaneous frequency (TDIF) [4], and (ii) different angular size of detection electrodes, including traditional wide (*e.g.* 90-120 degrees) and narrow (*e.g.*, as in NADEL ICR cells) configurations. Ion motion simulations performed using SIMION [2].

Results

Constant in time frequencies of ion motion, including cyclotron and magnetron components, have been

widely described in the FTMS literature using solutions of a well-known algebraic ICR frequency equation. However, the latter is applicable only when ions are trapped in a quadratic potential with axially-symmetric trajectories and their Coulombic interaction is negligible. However, in a real-life FTMS the instantaneous frequency of ion motion is normally time-dependent. Contrary to the quadratic trapping potential of the modern ICR cells, the non-quadratic trapping potential of the NADEL ICR cells appears to be responsible for the signal generation at the cyclotron frequency. By using the TDIF theory, in a combination with the experimental and computational results, we establish a relationship between the asymmetric trajectories of ion oscillations, transient modulation, detection electrode dimensions, and composition of resulting Fourier spectra. That allows us to formulate the analytical model for cyclotron frequency approach.

Conclusions

FTMS performance is limited by the space charge. It is responsible for peak interferences in mass spectra which result in peak shifts and, ultimately, coalescence. Experimental and simulation data with FT-ICR MS at the cyclotron frequency indicates it as a possible candidate for the space charge-tolerant FTMS. Presented here considerations and analytical models provide further insights into the fundamentals of this novel approach to FTMS.

Novel Aspect

An analytical model is developed for signal generation at the reduced and true cyclotron frequency in FT-ICR MS, including for the narrow aperture detection electrode ICR cells

References

1. Nagornov K. O., Kozhinov A. N., Tsybin Y. O., JASMS 28, 768–780 (2017)
2. Nagornov K. O., Kozhinov A. N., Tsybin Y. O., JASMS 29, 63-77 (2018)
3. Driver J. A., Nagornov K. O., Kozhinov A. N., Tsybin Y. O., Kharchenko A., Amster I. J., JASMS 29, 51-62 (2018)
4. Tsybin O. Y., Tsybin Y. O., Int. J. Mass Spectrom. 376, 75–84, (2015)

WP-189 / USE OF SIFT-MS IN THE ANALYSIS OF VOLATILE COMPOUNDS ARISING FROM THE ADDITION OF SODIUM HYPOCHLORITE AS AN IRRIGANT IN INFECTED ROOT CANAL TREATMENT

Kostas Ioannidis ⁽¹⁾ - **Sadia Niazi** ⁽¹⁾ - **Sanjukta Deb** ⁽¹⁾ - **Francesco Mannocci** ⁽¹⁾ - **David Smith** ⁽²⁾ - **Claire Turner** ⁽³⁾

King's College London, Dental Institute, London (1) - Transspectra Ltd, Newcastle-under-Lyme (2) - The Open University, School of Life, Health & Chemical Sciences, Milton Keynes (3)

Keywords: SIFT-MS, volatile organic compounds; dentistry

Introduction

Sodium hypochlorite (NaOCl) is a widely used disinfectant (1,2) and in dentistry, is routinely used as an irrigant in disinfecting infected root canals. In this study, we used selected ion flow tube mass spectrometry (SIFT-MS) to investigate the formation of toxic volatile organic compounds resulting from the reactions of NaOCl with a model system for infected root canals.

Methods

Samples representing a model tooth model system consisting of combinations of dentine powder (from ground freshly extracted molars), microbial suspensions and serum in a buffered suspension were incubated at 37°C under both aerobic and anaerobic conditions and then NaOCl was added. After incubation for 30 minutes, they were centrifuged. Volatile compounds in the headspace of the supernatants were analysed by selected ion flow tube mass spectrometry (SIFT-MS).

Results

During root canal chemomechanical preparation, NaOCl interacts with dentine, leading to denaturation of the organic matrix and changes in mineralised dentine (3,4). However, the production of potentially toxic organic compounds has not been widely investigated. SIFT-MS analysis of samples using the precursor ions H_3O^+ , NO^+ and O_2^+ showed that various volatile compounds were produced in abundance when NaOCl was added to different model tooth system components. NH_3 was produced when serum, dentine and bacteria were present, but there was much less in samples without serum. NaOCl also led to CH_3CN formation, possibly due to NaOCl's reaction with aldehydes and monochloramines, themselves formed through protein degradation by NaOCl. Furthermore, ions were detected at m/z 83, 85 and 87 using O_2^+ precursor. These ions are present with peak intensities in the ratio of 9:6:1, implying a 2 chlorine compound of which the only candidate is CHCl_2^+ . This closed shell ion is known to be the product of the O_2^+ reaction with chloroform (CHCl_3) (5).

Conclusions

SIFT-MS was an effective technique at identifying the volatile compounds produced in a model infected root canal system with the addition of the irrigant, NaOCl. Toxic compounds were formed by the interaction of NaOCl with the components of the model system, including chloroform and acetonitrile. These have the potential to contaminate the environment and lead to dental health worker exposure or may enter the patient's bloodstream.

Novel Aspect

The identification of toxic compounds produced during dental treatment is a novel application of SIFT-MS.

References

1. B. Wallace, M. Purcell, J. Furlong, J. Envir. Monit., 2002, 4, 35-42.
2. M. R. Khan, S. M. Wabaidur, Z. A. Alothman, R. Busquets, M. Naushad, Talanta, 2016, 152, 513-520.
3. M. Di Renzo, T. H. Ellis, E. Sacher, I. Stangel, Biomaterials, 2001, 22, 793-797.
4. A. Oyarzún, A. M. Cordero, M. Whittle, J. Endod., 2002, 28, 152-156.
5. P. Španěl P., D. Smith, Int. J. Mass Spectrom., 1999, 184, 175-181.

WP-190 / OPTIMAL SIGNAL-TO-NOISE AND PEAK ASSIGNMENT FOR TWO-DIMENSIONAL MASS SPECTROMETRY

Maria van Agthoven⁽¹⁾ - *Rémy Gavard*⁽¹⁾ - *Alice Lynch*⁽¹⁾ - *Tomos Morgan*⁽¹⁾ - *Christopher A. Wootton*⁽¹⁾ - *Lionel Chiron*⁽²⁾ - *Marc-André Delsuc*⁽³⁾ - *Mark P. Barrow*⁽¹⁾ - *Peter B. O'Connor*⁽¹⁾

University of Warwick, Department of Chemistry, Coventry (1) - CASC4DE, Research, Strasbourg (2) - IGBMC, Université de Strasbourg, Illkirch-Graffenstaden (3)

Keywords: Two-dimensional mass spectrometry, FT-ICR MS, ECD, ubiquitin, denoising

Introduction

Two-dimensional Fourier transform ion cyclotron mass spectrometry (2D MS) is a tandem mass spectrometry technique that requires no isolation and that maps the fragmentation patterns of all ion species in complex samples in a single 2D mass spectrum [1]. Scintillation noise in 2D mass spectra is minimized using the uncoiled random QR denoising algorithm (urQRd) [2]. The effect of denoising on the structural information available in 2D mass spectra is studied.

Methods

Seven 2D ECD mass spectra of ubiquitin were acquired on a 12 T solariX FT-ICR mass spectrometer. All spectra were processed using the SPIKE processing software with the urQRd denoising algorithm using

different ranks and without the urQRd denoising [3]. All horizontal fragment ion scans were peak-picked and the resulting peaklists were analyzed using the Themis software to assess the reproducibility of the 2D mass spectra [4].

Results

A 2D mass spectrum shows the fragmentation patterns of all the compounds in a sample. Each peak corresponds to a particular dissociation, where fragment ion scans, precursor ion scans and dissociation lines can be extracted from the 2D mass spectrum. The main source of noise is scintillation noise in the precursor ion dimension. The urQRd algorithm is routinely applied in the precursor ion dimension to lower the signal-to-noise ratio, but also removes peaks in the precursor ion dimension that are below the signal-to-noise ratio. In the fragment ion dimension, the signal-to-noise ratio is also drastically reduced, but many peaks are revealed by the orthogonal application of urQRd. We have extracted the peaklists from the fragment ion scans of seven 2D ECD mass spectra of ubiquitin processed with decreasing urQRd ranks, which corresponds to increasing levels of denoising. The peaklists have been compared using the Themis software to separate peaks corresponding to actual fragments from noise peaks.

Conclusions

This study enables a more confident assessment of signal-to-noise levels fragment ion scans in 2D mass spectra. With decreasing ranks of urQRd denoising, fragment ions are revealed in fragment ion scans. Below the optimal denoising level, fragment ions are removed alongside the noise and structural information is lost. This study improves the accuracy of the structural information that can be extracted from a 2D mass spectrum in terms of ubiquitin cleavage coverage.

Novel Aspect

Use of Themis software on protein top-down 2D MS. Study of the effect of urQRd denoising algorithm in 2D MS in the dimension orthogonal to the one it is applied to.

References

1. van Agthoven M. A., Delsuc M. -A., Bodenhausen G., Rolando C., Analytical and Bioanalytical Chemistry 405, 51 (2013)
2. Chiron L., van Agthoven M. A., Kieffer B., Rolando C., Delsuc M. -A., Proceedings of the National Academy of Sciences 111, 1385 (2014)
3. Chiron L., Coutouly M. -A., Starck J. -P., Rolando C., Delsuc M. -A., arXiv.org, e-Print Arch. Phys., 1-13 (2016).
4. Gavard R., Rossell D., Spencer S. E. F., Barrow M. P., Analytical Chemistry 89, 11383 (2017)

WP-191 / SIGNAL AMPLIFICATION MASS SPECTROMETRY BIOASSAY FOR ATTOMOLAR DETECTION OF GENES

Shuting Xu - Yu Bai - Huwei Liu

Peking University, College of chemistry, Beijing

Keywords: desorption electrospray ionization (DESI), mass tag, signal amplification, Wilms tumor gene

Introduction

The use of mass spectrometry (MS) as a readout system provides several advantages of high sensitivity, high resolution, rich structural information, and absolute quantification capability. Fast qualitative and quantitative analysis of biological samples, such as whole blood drops [1], tissue sections [2], and single cells [3], has been realized by MS, especially after ambient mass spectrometry (AMS) emerged [4,5]. However, it has limitations for ultrasensitive direct detecting of large biomolecules by AMS[6]. Herein, a

simple desorption electrospray ionization mass spectrometry (DESI-MS)-based bioassay utilizing signal amplification strategy was developed for ultrasensitive detection of large biomolecules, such as proteins and DNA.

Methods

A glass chip was functionalized with thiol-modified capture oligonucleotides, while gold nanoparticles were functionalized with both DNA capture oligonucleotides and mass tags. Here, a new ionic small molecule was synthesized as mass tag. The target, WT1 DNA (67 base long ssDNA oligonucleotide, an important biomarker for minimal residual disease (MRD) in leukemia)[7], was captured by a sandwich assay between biochips and nanoparticles. After being dried, the active biochips were detected by DESI-MS. The dissociated mass tags, together with internal standards, were detected by MS for quantitative analysis of target oligonucleotides.

Results

Interesting dissociation results were found that a large number of mass tags were dissociated from gold surface and produced abundant mass tag derivatives by DESI just in the open air without high-energy laser or high temperature. As hundreds of mass tags are carried by one gold nanoparticle, good amplification can be achieved for target molecule. Benefiting from the signal amplification, the detection of WT1 DNA could achieve ultrahigh sensitivity with a limit of detection of 0.5 amol (1 μ L, 0.5 pmol/L) and a wide linearity range from 1 amol to 500 amol.

Conclusions

A novel platform for large biomolecule detection was developed, which combined the signal amplification strategy and DESI-MS to achieve ultrahigh sensitivity and high efficiency detection of a disease biomarker (WT1 DNA). As a general method, this bioassay could be utilized for other DNA or protein biomarkers, which would be valuable for the early diagnosis in clinical diagnosis.

Novel Aspect

Multistep signal amplification with mass tags was realized by DESI-MS in the open air without high-energy laser or high temperature, and high sensitive detection (lower than attomole) of WT1 DNA was achieved by the simple DESI-MS bioassay.

References

1. Wang, H., Liu, J., Cooks, R. G., Ouyang, Z., *Angew. Chem. Int. Ed.* 49, 877-880 (2010).
2. Wiseman, J. M., Ifa, D. R., Song, Q., Cooks, R. G., *Angew. Chem. Int. Ed.* 45, 7188-7192 (2006).
3. Ferreira, C. R., Jarmusch, A. K., Pirro, V., Alfaro, C. M., Gonzalez-Serrano, A. F., Niemann, H., Wheeler, M. B., Rabel, R. A. C., Hallett, J. E., Houser, R., Kaufman, A., Cooks, R. G., *Reproduction Fertility and Development* 27, 621-637 (2015).
4. Takats, Z., Wiseman, J. M., Gologan, B., Cooks, R. G., *Science* 306, 471-473 (2004).
5. Cody, R. B., Laramée, J. A., Durst, H. D., *Anal. Chem.* 77, 2297-2302 (2005).
6. Zhang, C., Glaros, T., Manicke, N. E., *J. Am. Chem. Soc.* 139, 10996-10999 (2017).
7. Jorgensen, J. L., Chen, S. S., *Clin. Lymphoma, Myeloma Leuk.* 11, S49-53 (2011).

WP-192 / TWO-DIMENSIONAL (2D) FT-ICR MS AT HIGH RESOLUTION ON PRECURSOR IONS.

Christian Rolando⁽¹⁾ - ***Fabrice Bray***⁽¹⁾ - ***Julien Bouclon***⁽²⁾ - ***Philippe Pelupessy***⁽²⁾ - ***Marc-André Delsuc***⁽³⁾

CNRS, University of Lille, Villeneuve d'Ascq (1) - PSL, Ecole Normale Supérieure, Paris (2) - IGBMC, University of Strasbourg, Strasbourg (3)

Keywords: FT-ICR, 2 dimensional, high resolution on precursor ions

Introduction

Most of modern mass spectrometers use quadrupole analyzer isolation with a selection window in the range of one m/z to the best as the transmission is decreasing fast beyond. We developed also non-uniform sampling (NUS) acquisition for 2D FT-ICR MS which consists in skipping points randomly on the first dimension which give access to a higher resolution in the same amount of time.

Methods

Nevertheless, the maximum delay between the two encoding pulse increases as the resolution and reaches 1 second using for a 2 Megapoints acquisition with a 0.5 microsecond dwell corresponding on a 9.4 Tesla instrument to a resolution of 380,000 at m/z 200. Obtaining sounded 2D FT-ICR implies that ions stays in place without drifting during this time with a precision better than the ECD or IRMPD beam.

Results

We will describe here the optimization of the acquisition and of the treatment which allow us to obtain square 2D spectra. First of all we continue our analytical calculations of the ion movement to better understand the fundamental of 2D FT-ICR during the complete sequence. With this tool in hands the acquisition was optimized in two ways: (i) by modifying the Bruker IRMPD device to obtain a wider beam and a more precise control of the alignment; (ii) by a fine tune of the ICR cell for 2D acquisition (iii) and by using an optimized 2D pulse sequence. The data treatment was based on the iterative use of an improved version of the urQRd algorithm, which we developed previously for denoising, considering the missing points as noise.

Conclusions

Currently we are recording 2D of size F1 (precursor) 256k and size F2 (fragments) 256k points which allows us to reach a resolution of 50,000 at m/z 200 for precursors and fragments which is five times higher than our previous best results. The acquisition is performed with a NUS ratio 1/16 which means a physical size of 16k which can be recorded overnight.

Novel Aspect

We describe here for the first time the analysis of a complex mixture by MS/MS analysis without chromatography with high resolution on precursors.

References

1. Van Agthoven, M. A., et al. Towards analytically useful two-dimensional Fourier transform ion cyclotron resonance mass spectrometry. *Analytical and Bioanalytical Chemistry* 405, 51-61 (2013).
2. Bray, F., et al. Nonuniform Sampling Acquisition of Two-Dimensional Fourier Transform Ion Cyclotron Resonance Mass Spectrometry for Increased Mass Resolution of Tandem Mass Spectrometry Precursor Ions. *Analytical Chemistry* 89, 8589-8593 (2017).

WPS-S09 - Food additives and contaminants

Chairs: Jose Diana Di Mavungu, Roberta Galarini

WP-193 / DETERMINATION OF UNDESIRE SUBSTANCES IN EDIBLE OILS BY UHPLC-MS/MS

Alessandro Baldi Talini⁽¹⁾ - ***Yufeng Gao***⁽²⁾ - ***Zhenpeng Zhen***⁽²⁾ - ***Mingli Zhu***⁽³⁾ - ***Feng Qin***⁽⁴⁾ - ***Josh Ye***⁽⁴⁾

Perkinelmer In, Mass Spectrometry, Shelton (1) - Guangdong Provincial Bioengineering Institute (Guangzhou Sugarcane Industry Research Institute), Guangzhou, China (2) - Guangzhou Agricultural

Products Quantity and Safety Supervisory Institute, Guangzhou (3) - PerkinElmer Inc, Mass Spectrometry, Bolton (4)

Introduction

Edible oil production process may cause the formation of degradation products like polycyclic aromatic hydrocarbons (PAHs) or 3-MCPD esters and glycidol esters. Other contaminants include bisphenol A or phthalates from the packaging materials, mycotoxins, such as aflatoxins (AFs), ochratoxin A (OTA), deoxynivalenol (DON) and zearalenone (ZEA) and in some cases antioxidants which are frequently added to cooking oil for prolong shelf life. These undesired contaminants, if present in significant levels in oil, can pose serious health risk to consumers. Traditionally, these different classes of compounds are usually analyzed by different methodologies. Present study explores the possibility of analyzing all of the above compounds by a simple workflow of dSPE clean-up sample preparation followed by UHPLC-MS/MS analysis.

Methods

Oil sample obtained from local stores were extracted with methanol, followed by dispersive solid phase extraction (dSPE) clean-up. Analytes detection were via a highly sensitive PerkinElmer QSight 220 UHPLC-MS/MS system equipped with dual ESI/APCI sources.

Chromatographic separation was performed on a Brownlee™ SPP C18 (100x2.1 mm, 2.7 μm) column reversed-phase column using a gradient elution.

Under the optimized conditions, the method was validated in terms of its specificity, sensitivity, linearity, accuracy, and precision. Recoveries from sample matrix were evaluated by fortifications of analytes at three different concentrations.

Preliminary Results

Calibration curves were built by preparing several concentration levels of analytes standards in oil sample matrix (matrix-matched calibration) to overcome any matrix effects. The calibration curves obtained for all analytes were linear over three orders of magnitude with a R² value greater than 0.990. The LOQs for these analytes ranged from 0.1 to 500 μg/L.

The mean recoveries of those analytes from six samples spiked at different concentration levels were between 70% and 120%. The intra-day and inter-day variations, expressed as RSD, were less than 8%, respectively. The results demonstrated that the accuracy and precision of the present method were acceptable for routine monitoring purposes.

Identification and confirmation of the analytes in investigated real sample was achieved by comparing their respective retention times and peak area ratios between the quantifier and qualifier MRM transitions with those from the standards. Overall, the method developed in this study is simple, rapid, sensitive and robust. Demonstrated performance indicates that such unified method can be applied to determination of certain undesired substances in edible oil.

Future work would include further development of LC method for better separation which includes optimization of mobile phase conditions and MS parameters optimization for better sensitivity.

Novel Aspect

Simultaneous determination of different classes of undesired substances in edible oil by high pressure LC-MS/MS.

WP-194 / DEVELOPMENT AND VALIDATION OF AN LC-MS/MS METHOD TO QUANTIFY FUMONISINS AND HYDROLYSED FUMONISINS IN PORCINE SERUM

Gerlinde Bichl⁽¹⁾ - ***Tibor Czabany***⁽¹⁾ - ***Heidi-Elisabeth Schwartz-Zimmermann***⁽²⁾ - ***D. Schatzmayr***⁽¹⁾ - ***Ilse Dohnal***⁽¹⁾

Biomin, Research Center, Tulln (1) - University of Natural Resources and Life Sciences, Department of

Keywords: mycotoxin, fumonisin, serum, LC-MS/MS, ¹³C-labelled standards

Introduction

Fumonisin are mycotoxins that negatively impact animal and human health. The predominant form is fumonisin B1 (FB1), followed by FB2 and FB3. A novel feed additive, containing fumonisin esterase (FUMzyme®), can enzymatically remove the tricarballylic acid groups of fumonisins [1], dramatically reducing their toxicity [2]. We aim to quantify toxins and reaction products in biological matrices.

Methods

We describe an analytical method for the quantification of FB1, FB2, FB3 and all their fully (HFB) and partially (pHFBa/b) hydrolysed forms in pig serum. Due to their low bioavailability, the concentration of these twelve analytes in serum is in the low µg/L range. The sample extraction with acidified methanol-acetonitrile mix was optimised, prior to analysis by HPLC-ESI-MS/MS in MRM-mode.

Results

To compensate for losses during sample preparation and for matrix effects in ESI-MS/MS analysis, ¹³C-labelled internal standards (IS) proved necessary. They were generated by enzymatic hydrolysis of fully ¹³C-labelled FB1, FB2 and FB3.

Moreover, the chromatographic separation had to be fine-tuned, especially for pHFB2a/b and pHFB3a/b: FB2 and FB3 are positional isomers and fragment ions of sufficient intensity are not isoform-specific.

The optimised method was validated by spiking the twelve analytes into blank serum in ten different concentrations on three different days in triplicates.

The overall recoveries were between 85 and 110% for all analytes. In contrast, recoveries without IS-correction were only between 40 and 75% for FB2, FB3 and their pHFBs. The limit of quantification was between 0.14 and 0.67 ng/mL, except for FB1 and FB2 (2.9 and 1.5 ng/mL, respectively).

Conclusions

The quantitative determination of biomarkers in different biological matrices is a prerequisite for detailed analysis of animal health studies.

The validation data demonstrate the usability of the developed analytical procedure and measurement method for accurate identification and quantification of fumonisins and their hydrolysis products in pig serum.

Novel Aspect

A validated method to analyse FB1, FB2, FB3 and their hydrolysis products in serum, using enzymatically produced, ¹³C-labelled internal standards.

References

1. Heintl S., Hartinger D., Thamhesl M., Vekiru E., Krska R., Schatzmayr G., Moll W. D., Grabherr R., Journal of Biotechnology 145, 120 – 129 (2010).
2. Grenier B., Bracarense A. P. F. L., Schwartz H., Trumel C., Cossalter A. M., Schatzmayr G., Kolf-Clauw M., Moll W. D., Oswald I. P., Biochemical Pharmacology 83, 1465 – 1473 (2012).

WP-195 / IDENTIFICATION OF SYNTHETIC FOOD DYES BY NANO-ELECTROSPRAY IONIZATION

Lucie Borovcová - Eva Máčalová - Karel Lemr

Department of Analytical Chemistry, Faculty of Science, Palacký University, Olomouc

Keywords: dyes, food additives, mass spectrometry, nanoelectrospray ionization

Introduction

Food color is an important parameter making a product more attractive for customers. Food dyes, natural or synthetic, are commonly used for compensation of the loss of natural colors due to the processing or storage of food. Synthetic dyes are strictly regulated. In adulterated food amount of dyes can exceed permitted limits or dyes are not mentioned in the list of ingredients [1,2]. Identification of dyes by nanoESI applying simple sample preparation is presented.

Methods

Experiments were carried out using a Xevo TQD triple quadrupole mass spectrometer (Waters, Manchester, UK). The commercial inlet of instrument was replaced by a custom-made inlet with a sample cone and a heated capillary allowing nanoelectrospray ionization (nanoESI). Nanoelectrospray tips ($2 \pm 1 \mu\text{m}$ I. D., PicoTips emitter, New Objective, Woburn, USA) were filled with sample solutions. Standards of fifteen synthetic food dyes were obtained from Sigma Aldrich (Prague, Czech Republic).

Results

Food dyes were successfully ionized by nanoelectrospray ionization (nanoESI). NanoESI was applied for its tolerance to components of matrix and lesser contamination of a mass spectrometer due to significantly smaller flow rate of sample solution in comparison to electrospray. Mainly sodium adducts $[\text{M}+\text{Na}]^+$ or $[\text{M}-\text{Na}+2\text{H}]^+$ were observed in spectra of dyes in positive mode. Better ionization was observed in negative mode, typical ions $[\text{M}-\text{Na}]^-$, $[\text{M}-2\text{Na}+\text{H}]^-$ (generally $[\text{M}-n\text{Na}+(n-1)\text{H}]^-$) and $[\text{M}-3\text{Na}+\text{H}]^{2-}$, $[\text{M}-3\text{Na}]^{3-}$ were formed. Simple dissolution of dyes from real food samples was applied. For example, the presence of declared dye brilliant blue FCF was confirmed by nanoESI and MS/MS experiment in blue candy M&M's after its dissolution in methanol:water (1:1, v/v).

Conclusions

NanoESI allowed analysis of synthetic food dyes with simple and fast sample pretreatment. Typically better ionization efficiency was observed in negative mode. Applicability was confirmed by identification of dyes in real food samples.

Novel Aspect

NanoESI and MS/MS experiments represent a useful tool for simple identification of synthetic food dyes in food samples.

References

1. Gosetti F., Chiominatto U., Mazzucco E., Mastroianni R., Marengo E., Food Chemistry, 167, 454-462, (2015).
2. Huang H. Y., ChiuCh. W., SueS. L., Cheng Ch. F., J. Chromatogr. A, 995, 29-36, (2003).

Acknowledgment

The authors gratefully acknowledge the financial support from the Ministry of Education, Youth and Sports of the Czech Republic (LO1305) and Palacký University in Olomouc (IGA_PrF_2018_027).

WP-196 / ACRYLAMIDE DETERMINATION IN POTATO AND/OR CEREALS BASED FOOD BY LC-MS/MS METHOD: ONE YEAR OF MONITORING IN NORTH-WEST ITALY

Valentina Ciccotelli - Lucia Anna Masiello - Valentina Savio - Barbara Vivaldi - Carlo Ercolini

Istituto Zooprofilattico Sperimentale del Piemonte, Liguria e Valle D'Aosta, sezione di Genova, Genova

Keywords: acrylamide, LC-MS/MS, carbohydrate-rich food

Introduction

Acrylamide (2-propenamide) is a water-soluble molecule, which is formed by Maillard reaction during manufacturing or home-cooking of carbohydrate-rich food at high temperature [1], especially during browning process at temperature higher than 120°C, in food with high level of sugar and asparagine, like potato chips, French fries, crackers, bread, coffee and biscuit.

In 2002 the International Agency for Research on Cancer (IARC) classified Acrylamide in Group 2A as a probable human carcinogen [2], for this reason European Union established mitigation measures and reference levels for reducing the presence of acrylamide in food with the Regulation 2158/2017 [3].

The aim of this work was to develop and validate a method in LC-MS/MS coupled with QuEChERS approach to carry out a monitoring of processed foods in North-West Italy from 2017 to 2018.

Methods

An aliquot (2 g) of ground rehydrated sample was homogenized with 10 ml of acetonitrile. A prepared mix of QuEChERS was added and shaken for 10 min. The organic phase was separated by centrifugation. An aliquot of upper layer was dried in water bath at 45°C and 120 mbar. The sample was dissolved in a volume of Milli-Q water and filtered through a 0.22 µm PES membrane into vial, ready for analysis in LC-MS/MS.

Analysis was performed by an HPLC Accela 1250 pump (Thermo Scientific). Chromatographic separation was obtained by a reversed-phase Kinetex Biphenyl (100 x 2.1 mm ID, 2.6 µm) (Phenomenex). The eluents were 0.1% formic acid aqueous solution (eluent A) and 0.1% formic acid in methanol (eluent B), the elution was carried out in isocratic mode (total run 7 min), at flow-rate of 200 µl/min at 20°C. The injection volume was 5 µl. Mass spectral analysis was performed on a TSQ Vantage triple quadrupole mass spectrometer (Thermo Scientific), equipped with a heated ion spray interface (HESI). Two mass transitions (72.1/55.1 and 72.1/27.4) were monitored for acrylamide, providing unequivocal identification of it and only one mass transition (75.3/58.6) was used for the internal standard (acrylamide-d₃).

Method evaluation included: linearity, matrix effect, recovery, precision (as intra- and inter-laboratory RSD%) and LOQ, using certified reference materials (CRM) and spiked samples.

Results

Method validation was performed at three concentration levels: 25 µg/Kg, 120 µg/Kg and 625 µg/Kg. For the first level, we used a cereal based baby food spiked with a standard solution of acrylamide, we decided to fix the LOQ at this level; for the second level a biscuit based CRM and for the last level a potato based CRM, to cover all matrices of routine samples.

Linearity was evaluated over the concentration range from 25 µg/Kg to 1000 µg/Kg by constructing a calibration curve in solvent. Matrix effect was verified by comparing the calibration curve in solvent with calibration curve in matrix-post for the matrices: potato chips, biscuits, flakes and mashed potatoes.

Fourteen samples of different cereals or potato based matrices were extracted and analyzed. The occurrence of possible interferences was tested by plotting the acrylamide selected-ion chromatograms at the relative retention time interval ($\pm 2.5\%$) expected for its elution. Satisfactory selectivity is achieved when area peaks was under 10% of medium area peaks of sample at LOQ in the corresponding retention time windows for the quantification SRM transition.

Intra- and inter-day precision were determined by assaying six replicates of sample at three concentration levels in two different analytical batches. Repeatability was calculated by applying the Horwitz equation and the acceptable criteria was based on this equation: intra-repeatability RSD = 0.66*inter-repeatability RSD [3].

The recovery was evaluated on the same samples and ranged from 90% to 101%.

Fifteen samples of crackers, chips, baked potatoes, rice cream for children, cereals, biscuits, bread, rusks were purchased on the North-West Italian market. Baby food, bread and rusks showed values lower than LOQ; the highest concentration of acrylamide was found in potato chips (average concentration of 600 µg/Kg), the baked potatoes showed an average content of 269 µg/Kg; breakfast cereals, biscuits and crackers had concentrations from 27 µg/Kg to 54 µg/Kg.

Conclusions

The study describes a quick and easy method for acrylamide determination by QuEChERS approach coupled to LC-MS/MS. Validation of this method fulfilling recoveries (90%-101%), and showing a good linearity (from 25 µg/Kg to 1000 µg/Kg). Precision varies from 1.43% to 6.38% of intra-repeatability and from 3.73% to 10.36% of inter-repeatability.

The method was applied to food matrices such as crackers, chips, biscuits, bread, rusks and baby food. results confirm that acrylamide was found in ordinary consumer products, this underlines the importance of these monitoring to ensure consumer health.

Novel Aspect

There are relatively few methods in the scientific literature addressed to the analysis of acrylamide in processed food, baby food included. Our protocol distinctively combines simplicity of sample treatment but effective, applicability to very different matrices and compliance of the sensitivity requirements expected by European and Italian legislations.

References

1. Kim T. H., et al, Talanta, 131, 46-54 (2015)
2. De Paola E. L., Montevecchi G., Masino F., Garbini D., Barbanera M., Antonelli A., Food Chemistry, 217, 191-195 (2017)
3. Regulation (UE) n. 2158/2017 of 20 November 2017, Official J. of the Eur. Union L304, 24-44.

WP-197 / DETERMINATION OF PESTICIDES IN FRUIT, VEGETABLES AND FEED BY MULTIRESIDUE METHODS: ONE YEAR OF MONITORING IN NORTH-WEST ITALY

Valentina Ciccotelli⁽¹⁾ - *Raffaella Nappi*⁽²⁾ - *Barbara Rutigliano*⁽²⁾ - *Annamaria Galleggiante Crisafulli*⁽²⁾ - *Elisabetta Poma Genin*⁽²⁾ - *Enrico Alesso*⁽¹⁾ - *Barbara Vivaldi*⁽¹⁾

Istituto Zooprofilattico Sperimentale del Piemonte, Liguria e Valle D'Aosta, N. R. L. for pesticides in cereals and feed, sezione di Genova, Genova (1) - Istituto Zooprofilattico Sperimentale del Piemonte, Liguria e Valle D'Aosta, sezione di Cuneo, Cuneo (2)

Keywords: pesticides, cereals, fruit, vegetables, monitoring

Introduction

Pesticides are used in various combinations at different stage of cultivation or postharvest storage, to protect crops against a wide range of pests and fungi or to provide quality preservation. Numerous monitoring studies have demonstrated the presence of pesticide residue and their degradation products (metabolites) on crops, in environments and in the diet. These residue levels are usually low, below of effective action and interaction is not expected to occur. At the same time, literature data indicate that presence of pesticide mixtures in food and feed could be a potential risk to human and animal health.

The aim of this study was to monitor pesticides in fruit, vegetables and cereals through using of screening-quantitative multiresidue methods (GC MS/MS and UHPLC MS/MS) developed and validated for selected analytes. The survey was conducted in North West Italy (from 2016 to 2017) and these regions constitute one of the main cultivation areas of Italy. In these areas the main crops were cereal grains, rice, corn, fruit, vegetables and vineyards.

Methods

Screening quantitative multiresidue method (GC-MS/MS and UHPLC-MS/MS) were developed and validated for selected pesticides and their metabolites in fruit vegetable and cereals. The homogeneous sample was extracted with the help of acetonitrile. Samples with low content (<80%) require the addition of water before the initial extraction to get a total of approximately 10 g of water. After addition of

magnesium sulfate, sodium chloride and buffering citrate salts, the mixture was shaken intensively and centrifuged for phase separation. An aliquot of the organic phase was cleaned-up by dispersive solid phase extraction (D-SPE) employing bulk sorbents as well as magnesium sulfate for the removal of residual water. Following clean-up with amino-sorbents (*e.g.* primary secondary amine sorbent, PSA) extracts were acidified by adding a small amount of formic acid, to improve the storage stability of certain base-sensitive pesticides. The final extract was carried out by GC- and LC-based determinative analysis.

In addition of these method a multi-residue quantitative method for the analysis of 103 gas-chromatography amenable pesticides in dry cereal based animal feed was developed. The method entails a simple extraction of the rehydrated sample with acidified ethyl acetate. The organic phase was then separated by centrifugation, after the addition of anhydrous sodium sulphate, and filtered prior the final determination by GC-MS/MS.

Two mass transitions were monitored for each pesticide, providing unequivocal identification of contaminants. Where the method didn't include any clean-up step, the choice of the mass transitions to be monitored was a key-point of method development. Detected pesticides were then quantified using triphenyl phosphate as internal standard.

Method validation was carried out according to European guideline document SANTE 11945/2015 [1]. All the parameters studied (specificity, linearity, matrix effect, LOQ, recovery, precision and ruggedness) met the criteria set in the guideline at the three different spiking levels investigated (0,010 – 0,050 and 0,10 mg/kg). Also expanded uncertainty of the result was calculated using a top-down approach as indicated in the Codex Guidelines on estimation of uncertainty of results (CAC/GL 59-2006 amendment 2011).

Results

National monitoring results of the 2016 year have shown that nearly 48% of the foods sampled had no pesticide residues detected. Pesticide residues were detected in 52% of samples analyzed and the results showed a different distribution of active substances in the various commodity: 76% in fruit/vegetable, 33% in cereals and 80% in wine grapes/wine. The most frequently detected pesticide was fungicide and insecticide, followed by acaricide and herbicide. The study showed that pesticide residues exhibited a seasonal variation in terms of detection frequency and concentration, with spring and summer as critical seasons.

In particular N. R. L. for pesticides in cereals and feed analyzed 93 samples of feed and cereals, in 56 samples, from one to six pesticides were detected in quantifiable quantities (>LOQ 0.010 mg/Kg) but below the MRL (Maximum Residue Level).

The most common detected substances are: pirimiphos-methyl, most times associated with piperonyl-butoxide (used as preservative and not yet normated) and cis-deltamethrin.

While cypermethrin, chlorpyrifos-methyl, chlorpyrifos, endrin, malaoxon, desmediphamresmethrin, permethrin, azoxystrobin, fenazaquin and difeconazole were determined with less frequency. The most contaminated matrices are feed, rice and wheat.

Conclusions

Pesticide residues were detected in different matrices examined at values below to the legal limit. Active substances were identified in high percent in fruit/vegetable and wine/wine grapes with fungicide and insecticide mainly detected. However, for national application in monitoring programs it is necessary to develop new methods in multi and monoresidue screening which cover a broad spectrum of pesticides and metabolites.

References

1. Document N° SANTE/11945/2015 implemented by 01/01/2016

WP-198 / ROBUST MULTIRESIDUE PESTICIDE ANALYSIS IN BABY FOOD WITH ENHANCED SENSITIVITY GC-MS/MS TECHNOLOGY

Elena Ciceri⁽¹⁾ - ***Richard Law***⁽²⁾ - ***Tim Anderson***⁽³⁾ - ***Cristian I. Cojocariu***⁽²⁾

Thermo Fisher Scientific, CMD, Milan (1) - Thermo Fisher Scientific, CMD, Runcorn (2) - Thermo Fisher Scientific, CMD, Austun (3)

Keywords: Pesticides, food, GC-MSMS, TSQ9000

Introduction

Pesticides are widely used in food production and to control and manage various unwanted species. The detection, quantification, and confirmation of pesticide residue in food at trace levels require sensitive, selective and robust analytical instrumentation.

With ever increasing pressure to analyze contaminants at very low levels in a greater number of samples and with shorter turnaround times, laboratories seek continuous improvements in analytical instrumentation in order to maximize productivity and reduce downtime.

Methods

This study describes an optimized sample extraction and the analysis of multi-residue pesticides in baby food using femtogram-level sensitivity on a Thermo Scientific GC-MS/MS system. Sensitive detection and quantification of pesticides in baby foods is of particular importance as infants are more vulnerable to adverse health effects from these chemicals than adults. In this work, commercially available fruit and vegetable based baby-food samples were subjected to an optimized QuEChERS extraction method to isolate the pesticide residues using acetonitrile as the final extraction solvent. Direct analysis of extracts in acetonitrile is desired to avoid the need for solvent exchange to a more GC amenable solvent. The QuEChERS extracts were then analysed for pesticide residue content using a fast, targeted, timed-SRM method for >200 pesticides (inclusive of priority pesticides). The performance of the method focused on assessing various quantitative parameters for analysis of over 200 pesticides.

Results

The results show outstanding LODs (as low as 0.005 µg/Kg), robustness (over 900 matrix injections with no maintenance) and linear responses over up to 5 orders of magnitude (0.025 – 250 µg/kg) for the compounds investigated.

Conclusions

The data presented in this study demonstrate unprecedented method performance from sample preparation to sensitive and robust GC-MS/MS analysis in addition to automated data processing and reporting capabilities.

WP-199 / SENSITIVE AND ROBUST DETERMINATION OF PHTHALATES IN COOKING OIL USING ADVANCED ELECTRON IONIZATION SINGLE QUADRUPOLE GC-MS TECHNOLOGY

Elena Ciceri⁽¹⁾ - *Aaron Lamb*⁽²⁾ - *Cristian I. Cojocariu*⁽²⁾

Thermo Fisher Scientific, CMD, Milan (1) - Thermo Fisher Scientific, CMD, Runcorn (2)

Keywords: Pesticides, food, GC-MSMS, TSQ9000

Introduction

Phthalates are emerging food contaminants linked to endocrine disruption through ingestion of foods exposed to these common plasticizers. Due to their endocrine mimicking properties, phthalates are considered as priority pollutants and various countries around the world have already introduced maximum acceptable levels in food contact materials (FCMs) for compounds such as DMP, DEP, DAP, DIBP, DBP, DXHP, BBP, DCHP, DEHP, DINP, DNOP and DIDP.

Phthalates are lipophilic in nature and the propensity for these compounds to leach from FCMs into food

products such as cooking oils is significant. Oils are complex mixtures of triacyl glycerides that are difficult to chromatograph and present a challenge in terms of selectivity, sensitivity and robustness for GC-MS analysis of phthalates.

Methods

With the innovative design of the new advanced electron ionization source however less frequency of source cleaning is required as the axial ionization in this source leads to increased ionization efficiency. This means the source filament can be operated at a reduced emission current which in turn means less ionization of complex matrices in the source reducing the cleaning interval.

Results

The enhanced sensitivity of the new source down to low ppt levels (0.01 ng/mL) also means that the sample matrix can be diluted more or the split ratio can be increased which further reduces the amount of contamination in the GC flow path.

Here the use of a novel AEI high performance electron ionization (EI) source in conjunction with a new single quadrupole mass spectrometer provides the sensitivity and selectivity to analyze phthalates in cooking oil at low part per trillion (ppt) levels with low liquid injection volumes.

Conclusions

In this work the Thermo Scientific™ ISQ 7000™ GC with Advanced Electron Ionization Source mass spectrometer was evaluated for phthalate selectivity, sensitivity, recovery, repeatability, robustness and linearity of response.

WP-200 / BEE POLLEN AND PYRROLIZIDINE ALKALOIDS: A CONTAMINATION PROFILE ANALYSIS USING LC-MS/MS

Luciana De Jesus Inacio - Lorena Lucatello - Vittoria Bisutti - Ambra Sartorato - Roberta Merlanti - Francesca Capolongo

University of Padua, Department of Comparative Biomedicine and Food Science, Viale dell'Università 16, 35020 Legnaro, (Pd) Italy

Keywords: Pyrrolizidine alkaloids, bee pollen, LC-MS/MS, food safety

Introduction

Pyrrolizidine alkaloids (PAs) and their N-oxides (PANOs) contamination in foodstuffs and dietary supplements, such as honey and bee pollen, represents a potential health risk for consumers. This study aimed the development and validation of an LC-MS/MS method to quantify PAs/PANOs levels in bee pollen samples coming from local beekeepers (LB; Veneto/ Italy), great retail and web (GR).

Methods

The developed method was validated according to the Commission Decision 2002/657/EC. For this purpose, blank bee pollen samples were spiked with 17 PAs/PANOs standard solutions at three levels (2, 10 40 µg/Kg). PAs and PANOs were extracted from pollen using PLEXA cartridges (200 mg/6 mL, Agilent, USA) and the LC-MS/MS analysis was conducted on an LTQ-XL Ion trap operating in positive ion mode.

Results

All results obtained for recovery, linearity, repeatability and reproducibility in the validation study were in agreement with the limits reported by Decision 2002/657/EC. The limit of quantification was 0.4 µg/kg for all analytes. The validated method was applied to 81 bee pollen samples, of which 50 were positive: 20 were from LB and 30 from GR. The average total concentration of PAs/PANOs in samples from LB and GR

were 27.8 µg/kg and 470.9 µg/kg, respectively. The identified PAs/PANOs could be categorized into three classes and their average concentrations were as follows: Eupatorium-Type 42.90 µg/kg (LB) vs. 162.76 µg/kg (GR); Echium-Type 7.46 µg/kg (LB) vs. 465.54 µg/kg (GR); Senecio-Type 16.24 µg/kg (LB) vs. 176.62 µg/kg (GR). Most samples (17) from LB presented Senecio-Type PAs/PANOs, whereas most samples (24) from GR presented Eupatorium-Type PAs/PANOs.

Conclusions

The validated LC-MS/MS method showed to be sensitive and suitable for the quantification of the 17 selected PAs/PANOs markers in bee pollen. Samples from LB and GR presented a very different PAs/PANOs contamination profile reflecting their different botanical and geographical origin. The average total concentration of PAs/PANOs was 17 times lower in samples from LB than in those from GR.

Novel Aspect

This study should highlight any hazard due to the import/export of pollen and suggest management strategies to beekeepers to ensure a healthy product.

References

1. European Commission Decision 2002/657/EC, Off J Eur Communities, L 221/8:8–36 (2002).
2. Kempf M. et al., Mol Nutr Food Res. 54, 292–300 (2010).
3. Dübecke A., Beckh G., Lüllmann C., Food Addit Contam Part A, 28 (3), 348–358 (2011).
4. Mulder P. P. J. et al., EFSA Supporting Publ. EN-859, 1–114 (2015).
5. Kast C. et al., Food Addit Contam Part A, 35 (2), 316–327 (2018).

WP-201 / TOWARDS UNTARGETED CHEMICAL FOOD SAFETY ASSESSMENT

Grégoire Delaporte - Mathieu Cladière - Valérie Camel

UMR Ingénierie Procédés Aliments, AgroParisTech, Inra, Université Paris-Saclay, 91300 Massy, France

Keywords: UPLC-HRMS; foodomics; contaminants; chemometrics; trace analysis

Introduction

Ensuring food safety remains a major public health issue. Recent food safety crisis (*i.e.* eggs contaminated by fipronil in Europe) clearly illustrated the limits of analytical targeted approaches. The need for untargeted approaches has recently merged, especially for detecting non expected contaminants [1]. These new approaches raise many analytical challenges, due to contaminants diversity, low-levels concentrations (µg/kg) and food matrices complexity.

Methods

An untargeted foodomics approach was developed for blind food chemical contaminants detection. It relies on four pillars from metabolomics and chemometrics: (i) a broad-range UHPLC-HRMS method based on minimal sample treatment; (ii) feature detection based on XCMS and data cleaning tools; (iii) blind discrimination of sample groups by Independent Component Analysis; (iv) powerful data exploration strategies for annotation of suspect compounds.

Results

Tea has been chosen as a development food matrix since it is the most consumed hot beverage over the world, offering a high diversity and complexity as well. The proof-of-concept was initially developed on a rather simple case (known contaminants characterizing the potential diversity at different levels with a homogeneous matrix), and then successfully applied on more complex cases, such as the detection of an unknown contamination or the simultaneous analysis of samples from different brands of black and green

tea. Efficient data filtration, normalization and scaling strategies are key issues since a discerning choice of tools enables a blind discrimination between control and spiked sample groups at levels as low as 10 µg/kg, as well as the detection of unexpected contaminants. Annotation step appears to be the main challenge, but in-house data mining tools are key factors to reduce data complexity and also provide the user a quick feedback about potential identity or property (presence of an isotope or potential adduct) of unexpected compounds.

Conclusions

In this work, a metabolomics-like strategy was set up to detect and characterized a contamination in complex food samples at level down to 10 µg/kg. On top of that, the approach, based on open and free tools, was shown to have the capability to detect unexpected contaminants, which is critical for future chemical food safety assessment. No doubts that our results will encourage the development of such strategies for chemical food safety assessment within the next years.

Novel Aspect

We propose an original methodology for untargeted food safety assessment (down to 10 µg/kg) validated on several contamination scenarios for a complex real food matrix (tea).

References

1. J. P. Antignac, F. Courant, G. Pinel, E. Bichon, F. Monteau, C. Elliott, B. Le Bizec, *TrAC - Trends Anal. Chem.* 30 (2011) 292–301. doi:10.1016/j.trac.2010.11.003.

WP-202 / CHARACTERISATION AND QUANTIFICATION OF KNOWN/UNKNOWN MIGRATING SUBSTANCES IN FOOD USING LC & GC HIGH RESOLUTION MASS SPECTROMETRY

Yves-Alexis Hammel - Baptiste Bodet - Sander Koster

Nestle Research Center, Institute Food Safety & Analytical Science, Lausanne

Keywords: Food Contact Material, Mass Spectrometry, semi-quantitation, food simulant

Introduction

A method was developed that enables to 1) detect and categorise known and unknown migrating substances from packaging into food, 2) to semi-quantify these known and unknown substances using CAD and FID and 3) define an exposure threshold of unknown substances above which these will be identified. The methodology is applied to packaged food in comparison with measurements performed in so-called food simulants.

Methods

The developed approach consist of a QuEChERS extraction applied to the food sample and a migration with food simulant of the corresponding empty packaging. Then the extracts are analysed using GC-FID/PCI-HRMS and LC-CAD-HRMS. A semi-quantitation based on FID and CAD detection is applied. LC-HRMS and GC-PCI are used for identification purpose.

Results

A process was developed that enables to classify detected substances into known and unknown substances that migrate from the packaging material into food. Using our quantification/exposure based method, it is possible to focus on substances that are relevant to be considered for risk assessment.

Conclusions

The strategy for sample preparation of unknown substances will be presented together with measurements in food and food simulant. Levels of migrating substances in food will be compared with those in simulant

which helps to mitigate the concern of these migrating molecules in the risk assessment of packaging materials.

Novel Aspect

The approach allows to analyse known/unknown migrating substances in food giving real levels of migrating molecules in contrast to approaches used today that rely on food simulants.

WP-203 / MERCAPTOMIC: TOWARDS THE PROFILING OF REACTIVE METABOLITES IN TOXICOLOGY

Emilien Jamin⁽¹⁾ - *Robin Costantino*⁽¹⁾ - *Jean-François Martin*⁽¹⁾ - *Laurent Debrauwer*⁽¹⁾ - *Françoise Guéraud*⁽²⁾

INRA, MetaboHUB-MetaToul-AXIOM platform, Toulouse (1) - INRA, Toxalim, Toulouse (2)

Keywords: Metabolomics, Untargeted Mercapturic acid, 4-HNE

Introduction

The link between food-toxic effect and disease occurrence requires, the most complete characterization of exposure. In food safety, consumption data crossed with food contamination data or biomonitoring data allow evaluating exposure only in a targeted way. Based on our previous results in exposomic [1], food or environmental toxicology should focus on mixture of compounds, mostly at low doses, and in an untargeted way to identify unknown compounds.

Methods

We developed an untargeted method to profile electrophilic metabolites which are most of the time toxic, due to their chemical reactivity towards bio-molecules such as DNA. The main line of cellular defense against that is the conjugation to the tripeptide glutathione, which is metabolized in mercapturic acid conjugate to be excreted into urine. Therefore, mercapturate conjugates were monitored through their characteristic neutral loss observed in MS/MS.

Results

This approach was applied to the study of different groups of rats fed diets containing various oils. According to our previous results on lipid peroxidation [2] these diets led to the production of different aldehydes conjugated to mercapturic acid. The most famous is DHN-MA which corresponds to the mercapturate conjugate of 4-hydroxynonenal (4-HNE), which is commonly used as a biomarker of lipid peroxidation [2]. Using our methodology, we were able to detect without a priori, dozens of mercapturate conjugates, including DHN-MA and other known conjugated aldehydes. Furthermore, our approach also allowed the detection of conjugates of unexpected aldehydes, and of other chemical classes, for which putative identifications have been proposed based on complementary structural analyses. Interestingly, multivariate statistical analyses of the HRMS signals carried out on the mercapturate conjugates yield a better characterization of the studied animal groups compared to results obtained from a classic untargeted metabolomic approach.

Conclusions

Using the “all ions MS/MS” mode (MSE) of a Synapt G2-Si mass spectrometer it is possible to detect all the metabolites displaying the characteristic loss of a mercapturate conjugate, which represents a reactive metabolites. Thus, among these numerous metabolites representative of our exposure, it seems possible by our approach to focus on toxic ones.

Novel Aspect

We developed an untargeted profiling method of mercapturate conjugates, which are representative of

reactive metabolites in our exposome.

References

1. Jamin E. L. et al., Anal Bioanal Chem, 406, 114-1161 (2014).
2. Guéraud F., Free Radic Biol Med, 111, 196-208 (2017).

For information please contact: emilien.jamin@inra.fr

WP-204 / FAST AND AUTOMATIC DATA PROCESSING FOR TARGETED AND NON-TARGETED PESTICIDE ANALYSIS

Jasmin Meltretter⁽¹⁾ - ***Jianru Stahl-Zeng***⁽¹⁾ - ***Vanaja Raguvaran***⁽²⁾ - ***David M. Cox***⁽²⁾ - ***Julia Jasak***⁽¹⁾

SCIEX, Darmstadt, Germany (1) - SCIEX, Concord, Canada (2)

Pesticides are widely used in agriculture, growing of fruit and vegetables, and in lots of garden plots. The main aims are to protect the crops from pests and to control the growth of unwanted plants, both to minimize crop loss. The use of pesticides is strongly regulated to ensure that any hazards for humans, animals, and the environment are minimized. To monitor pesticide residues in food, LC-MS/MS in combination with QuEChERS sample preparation is today the method of choice for the majority of the relevant pesticides. With the SANTE/11813/2017 document, the European Commission published specific guidelines on quality control and validation procedures for analytical methods for pesticides in food and feed. In order to analyze more than 1000 pesticides, powerful data processing software, which ensures fast and easy data review, is strongly demanded by routine food testing labs.

The new SCIEX OS 1.4 software makes the data processing fast and easy. It includes powerful features like the automatic outlier removal feature, which will automatically fit the regression type and remove standard outliers from the calibration curve based on regression fit, accuracy, precision of replicates, and a user-defined outlier tolerance. This automatism tremendously reduces the time to check the calibration curves for all analytes.

Furthermore, the pesticide data have to be reviewed in accordance to the SANTE/11813/2017 document. All identification requirements of this guideline for both unit mass resolution and accurate mass measurement including retention time, mass accuracy of two ions and their specific ion ratio can be easily and fast reviewed in the software using a user-friendly traffic light system.

The SCIEX OS 1.4 software supports both QTOF data acquired with the SCIEX X500R QTOF system and the SCIEX TripleTOF® system and unit mass resolution data acquired with a SCIEX QTRAP® or SCIEX Triple Quad™ system.

WP-205 / SIMULTANEOUS DETERMINATION OF PROPOXYCARBAZONE AND ITS METABOLITE IN FOOD BY UHPLC-MS/MS, APPLYING A QUECHERS EXTRACTION BASED METHOD

Marta Vargas Pérez - ***Jesús Marín Sáez*** - ***María Luisa Jiménez Medina*** - ***Francisco Javier Egea González*** - ***Antonia Garrido Frenich***

Research Group "Analytical Chemistry of Contaminants", Department of Chemistry and Physics, Research Centre for Agricultural and Food Biotechnology (BITAL), University of Almeria, Agrifood Campus of International Excellence, ceiA3, E-04120 Almeria, Spain.

Keywords: Propoxycarbazone, metabolite, 2-hydroxy-propoxycarbazone, QuEChERS, UHPLC-MS/MS

Introduction

Nowadays, it is of great importance the monitoring of parent compounds as well as potential metabolites,

otherwise the effects of the use of chemicals would not be fully studied. The analysis of pesticide residues including metabolites in food matrices is a challenge, especially due to the low concentration of analytes and the large quantities of interfering substances that can be coextracted.

Methods

For the extraction of the compounds, several modifications were performed to the acetate-buffering QuEChERS method [1] depending on the commodities tested (lettuce, beetroot, honey and soybean meal). For analyte determination, ultra high performance liquid chromatography coupled to a triple quadrupole mass spectrometry analyzer was used. In addition, the cleaning stage was optimized, evaluating the effect of several sorbents, as well as the effect of the matrix during the quantification stage.

Results

The suitability of the method for the analysis of the target compounds in the different commodities was properly validated according to the SANTE document. [2]Linearity was evaluated using matrix matched standard calibration by analyzing spiked blank samples of each of the matrices under study at six concentration levels (from 5 to 200 µg/kg). R² was 0.99 in all the tested matrices. Lettuce and beetroot did not show significant matrix effect; however, honey and soybean meal have shown a strong matrix-effect. Recovery and precision was studied at two concentration levels (10, 50 µg/kg), obtaining adequate recoveries (73-110%), intra-day (4 - 19%) and inter-day precision (5 - 20%). The limits of detection (3-6 µg/kg) and quantification (10-25 µg/kg) were lower than the maximum residue limits (MRLs)[3] in the matrices evaluated. The validated method was applied to the analysis of 20 samples, but residues of target compounds were not detected in any of the samples analyzed.

Conclusions

The modified QuEChERS approach followed by clean-up in combination with liquid chromatography coupled to triple quadrupole method provided suitable selectivity and sensitivity for the simultaneous determination of propoxycarbazone and its metabolite in lettuce, beetroot, soybean meal and honey. The method is appropriate for routine monitoring of residues in vegetables, cereals and food samples to ensure food safety.

Novel Aspect

This is the first method that includes the simultaneous determination of propoxycarbazone and its metabolite in food stuffs, as it is established by regulations.

References

1. AOAC Official Method 2007.01, Pesticide residues in foods by acetonitrile extraction and partitioning with magnesium sulfate, 2007.
2. Guidance document on analytical quality control and method validation procedures for pesticide residues and analysis in food and feed. SANTE/11813/2017. European Commission. Directorate General for Health and Food Safety.
3. Peer review of the pesticide risk assessment of the active substance propoxycarbazone (variant evaluated propoxycarbazone-sodium). The EFSA Journal, 14, 4612-4637, (2016)

WP-206 / SENSITIVE DETERMINATION OF 44 TOXIC PYRROLIZIDINE ALKALOIDS IN SPICES AND CULINARY HERBS USING LIQUID CHROMATOGRAPHY TANDEM-MASS SPECTROMETRY

Florian Kaltner⁽¹⁾ - ***Michael Rychlik***⁽²⁾ - ***Manfred Gareis***⁽¹⁾ - ***Christoph Gottschalk***⁽¹⁾

Chair of Food Safety, Ludwig Maximilian University of Munich, Oberschleissheim (1) - Chair of Analytical Food Chemistry, Technical University of Munich, Freising (2)

Keywords: pyrrolizidine alkaloids, spices, culinary herbs, method development, tandem mass spectrometry

Introduction

1,2-dehydropyrrolizidine alkaloids (PA) and their *N*-oxides (PANO) are phytotoxins with hepatotoxic and mutagenic potential. PA presumably occur in plant species worldwide and are responsible for various recorded food and feed poisonings. Main source of consumers' exposition is food contaminated with seeds, pollen, nectar or parts of PA plants. Especially (herbal) tea infusions, food supplements, honey and, recently, spices can be contaminated with PA [1].

Methods

In total, 207 samples of spices and culinary herbs from Europe, Asia, Africa and America were investigated for occurrence of overall 44 PA. Therefore, a sensitive LC-MS/MS method was developed and validation was done for each group of spice and culinary herb matrices. Samples were extracted under acidic conditions and extracts were cleaned using strong cation exchange cartridges. Quantification was performed via external matrix-matched calibration.

Results

The method allowed the analysis of 44 PA within a total runtime of 16.0 min after a 6.5 min equilibration of the column. Limits of detection (LOD) and quantification (LOQ) ranged from 0.01 to 0.17 µg/kg and from 0.03 to 0.54 µg/kg, respectively. Precision of the method, expressed as the coefficient of variation of a six-fold determination of respective samples, was 0.9 % and 8.0 % for all analytes, except merepoxine and merenskinine and their *N*-oxides. Obtained recoveries lay between 45.1 % and 118.4 %, only single PA showed outlier data.

After the investigation of 151 of the 207 samples, in 60.3 % at least one PA was found. The obtained sum levels of PA-contamination ranged from 0.13 µg/kg in a rosemary sample to 24.62 mg/kg in an oregano sample. Including all positive samples, mean PA-sum level was 967.62 µg/kg. The most abundant detected PA in the samples were europine/*N*-oxide, heliotrine/*N*-oxide and lasiocarpine/*N*-oxide, which occur in plants of *Heliotropium* spp., a genus of the Boraginaceae family [2].

The final results will be reported.

Conclusions

A LC-MS/MS method for the analysis of 44 PA in spices and culinary herbs was developed and validated. The applicability was successfully proven by investigating 207 samples of different origins. The study showed that spices could contribute to the overall exposition to PA/PANO. However, this contribution is generally quite low due to the small daily intake of herbs and spices in comparison to other food commodities.

Novel Aspect

The developed method is the first target LC-MS/MS method being able to quantify 44 different PA and which was applied to a representative number of spice and culinary herb samples.

References

1. Federal Institute for Risk Assessment, Opinion No.030/2016, Berlin (2016).
2. Boppré M., Food Additives & Contaminants: Part A, 28, 260-281 (2011)

WP-207 / ANALYTICAL METHOD DEVELOPMENT OF MULTI RESIDUE PESTICIDES IN LIVESTOCK PRODUCTS BY GAS CHROMATOGRAPHY-MASS SPECTROMETRY(GC-MS/MS)

Hye Soon Kang - Eun Jeong Kim - Moo-Song Lim - Eun Jung Chang - KiMi No - Dong me Choi - Hyun-cheol Choi

Division of Hazardous Substances Analysis, Gyeongin Regional Office of Food and Drug Safety, 137 Juanro, Nam-gu, Incheon, 22133, Republic of Korea

Introduction

The Gas chromatography-mass spectrometry (GC-MS/MS) multi-residue method for the analysis of 55 pesticides in livestock products has been developed. The monitoring of pesticide residues in food is nowadays a priority objective in order to get extensive evaluation of food quality and to avoid possible risks to human health. The pesticides included in this study belong to different chemical families of insecticides, herbicides, fungicides and so on. GC-MS is rapidly becoming an accepted technique in pesticide residue analysis for regulatory monitoring purposes. This survey was aimed to establish for determination of multi-residue method for 55 pesticides in commercial livestock products by GC-MS/MS.

Methods

Briefly, the homogenized samples were mixed with 1% acetic acid in acetonitrile. The mixture was vigorously shaken, then added sodium chloride and $MgSO_4$ repeated the shaking process, and cleanup. After a centrifugation step, the resulting extract with solid-phase extraction using Florisil cartridge. The extracts were transferred to vials for analysis by GC-MS/MS. Pesticides were confirmed by their retention time and multiple reaction monitoring of three fragment ions by GC-MS/MS.

Results & Conclusions

Recovery, precision, accuracy, linearity, limit of determination (LOD), and limit of quantification (LOQ) in the analytical method were validated in different matrices (Pork, Beef, Chicken, Egg, Milk). Calibration curves were linear for most of the compounds studied. The recoveries obtained at fortified levels of 0.001~0.05 mg/kg were 60~120%, with relative standard deviations (RSDs) of $\leq 30\%$, with a few exceptions. This method we have developed can quickly and accurately analyze 55 kinds of pesticides used in livestock products. The methods will be proposed in the Korean Food Standards Codex as the analysis of residual pesticides in livestock products.

Novel Aspect

A GC-MS/MS multi-residue method for the simultaneous analysis of 55 kinds pesticides in livestock products.

WP-208 / LIMITATIONS OF GAS CHROMATOGRAPHY–SINGLE QUADROPOLE MASS SPECTROMETRY IN THE IDENTIFICATION OF COMPOUNDS MIGRATING FROM FOOD CONTACT MATERIALS

Csaba Kirckeszner - Zsuzsanna Eke - Zoltán Nyiri - Noémi Petrovics

Joint Research and Training Laboratory on Separation Techniques, Eötvös Loránd University, Budapest

Keywords: food contact materials, migration, GC-EI-Q-MS, GC-EI-TOF-MS, non-intentionally added substances

Introduction

Plastics are commonly used for food packaging. From these plastics potentially hazardous compounds can migrate into the food. It is difficult to establish the chemical structure of migrated compounds especially for the non-intentionally added ones. They are often detected with GC-EI-Q-MS and identified using NIST database without any further confirmation [1]. But to what extent can such results be confirmed?

Methods

The ability of GC-EI-Q-MS to identify migrating compounds was tested using plastic food contact materials (FCM). Degradation and migration into a food simulant was assisted with microwave radiation. The extracts

were analyzed by GC-EI-Q-MS. Preliminary identification of the migrants was based on the obtained mass spectra and a NIST library search. Confirmation was attempted using GC-EI-TOF-MS measurements and analytical reference materials, if available.

Results

As a result of the microwave irradiation of the polypropylene samples plenty of compounds were formed and migrated to 90 V/V% isooctane and 10 V/V% ethanol (solution for simulating fatty food). With GC-EI-Q-MS at least four compounds could be detected for all 26 samples. Using the NIST mass spectral database 11 different compounds could be preliminary identified with good matching factors (>80%). Out of these, 8 could be confirmed with GC-EI-TOF-MS with less than 5 ppm mass difference for all occurrences. On the other hand the Irgafos 168 – a commonly used antioxidant - and a degradation of it was found in 24 cases with relatively low matching factors (40-55%). All these occurrences could however be confirmed with both GC-EI-TOF-MS and analytical reference materials. Furthermore another common antioxidant, the butylated hydroxyl toluene (BHT) was first identified in 12 samples, but always under 70%. All these samples also contained its degradation product that was consequently misidentified with the NIST library search.

Conclusions

Our results prove that the GC-EI-Q-MS is a powerful tool in the detection of compounds migrating from FCM. However, identification of these compounds cannot reliably be based solely on a NIST mass spectral database search even though it is a rather widely used approach. Reliability of the identification can notably be improved with GC-EI-TOF-MS and analytical reference materials.

Novel Aspect

The identification of semi-volatiles migrating from FCM is often based on mass spectra of a single type of MS. Reliability of this approach is questionable.

References

1. C. Nerin, P. Alfaro, M. Aznar, C. Domeno, *Analytica Chimica Acta*, 775, 14–24 (2013).

WP-209 / MULTI-ANALYTE METHOD BASED ON LIQUID CHROMATOGRAPHY COUPLED WITH SINGLE QUADRUPOLE MASS ANALYZER FOR FAST AND EFFECTIVE CONTROL OF FOOD ADDITIVES

Ales Krmela - Vera Schulzova - Jana Pulkrabova - Jana Hajslova

UCT Prague, Department of Food analysis and nutrition, Prague

Keywords: food additives; liquid chromatography; mass spectrometry

Introduction

Food additives are widely employed by food industry to provide number of benefits. The use of food additives is regulated, for some of them maximum limits have been established. To control safe use of food additives and protect consumers' health, effective methods are needed. In line with current trends, this study was performed to replace number of 'conventional methods' based on LC-UV by multi-analyte method employing MS detector with single quadrupole analyzer.

Methods

Presented method uses an ultra performance liquid chromatography (UPLC Acquity, Waters) coupled to a simple single quadrupole MS (QDa, Waters) using a gradient of 5 mM aqueous solution of ammonium acetate and methanol while separating analytes on reversed phase BEH C18 column (Waters).

Results

The developed method enables simultaneous determination of number of synthetic food additives fairly

differing in physicochemical properties. Synthetic dyes, sweeteners, preservatives and several other compounds, vanillin, ethyl vanillin, theobromine and caffeine, were involved on the analytes list. Method was validated for sweet drinks and fruit concentrates. Repeatability, expressed as relative standard deviation (RSD), is less than 10 % for all analytes with recovery ranging from 80 to 120 %. Limits of quantification range from 0.1 to 1 µg/mL which is still considerably lower than legislation limits for food additives which are often ten or hundred times higher. Considering all, QDa seems to be a suitable detector for this kind of analysis. Method was applied for the analysis of fruit concentrates and fruit wines revealing a wide use of food additives.

Conclusions

QDa MS is an appropriate tool for rapid screening of synthetic food additives of which analyzed samples contained a significant amount despite not exceeding legislative limits.

Novel Aspect

Presented method enables simultaneous determination of 19 synthetic food additives fairly using a relatively cheap and simple single quadrupole MS.

WP-210 / AFLATOXIN M1 RESIDUES DETECTION IN MILK, INFANT FOOD AND CHEESE BY LC-MS/MS

Gian Luca Ferro⁽¹⁾ - *Stefania Massafra*⁽¹⁾ - *Paola Stella*⁽¹⁾ - *Enrico Pennazio*⁽¹⁾ - *Marta Leporati*⁽²⁾ - *Daniela Marchis*⁽²⁾ - *Marilena Gili*⁽¹⁾

Istituto Zooprofilattico Sperimentale del Piemonte, Liguria e Valle d'Aosta, S. S. Ricerca Residui, Torino (1) - Istituto Zooprofilattico Sperimentale del Piemonte, Liguria e Valle d'Aosta, S. S. Igiene dei Mangimi, Torino (2)

Keywords: Aflatoxin M1, milk, food safety, LC-MS/MS

Introduction

Aflatoxins are mycotoxins produced by some species of *Aspergillus*. In dairy cows, ingested Aflatoxin B1 is metabolized into Aflatoxin M1 (AFM1), which is eliminated through milk. Due to its high hepato carcinogenic potential, the level of AFM1 permitted is strictly regulated: EU Commission has established a maximum limit of 0.050 µg/kg for AFM1 in milk and a more restrictive limit of 0.025 µg/kg in milk and food intended for consumption by nursing infants and children[1].

Methods

First, a liquid-liquid extraction (LLE) was performed on 20 g sample. As extraction solvent, ethyl acetate was chosen for milk and dichloromethane for infant food, dessert and cheese. A delipidation was needed for creamy infant food and cheese. LC-ESI-MS/MS analysis was carried out using a Zorbax SB Phenyl RRHT column (2.1 x 50 mm, 1.8 µm). The method was fully validated as screening assay; the validation as confirmatory quantitative analysis is in progress.

Results

Milk extraction was performed as reported by Biancardi et al [2]; regarding to infant food, dessert and cheese, LLE has to be performed by homogenization with 100 ml dichloromethane and 10 g diatomaceous earth; after filtration, 20 ml of raw extract was dried under a stream of nitrogen; the residue was dissolved in HCOOH 0.1%:CH₃CN (90:10, v/v) and defatted with hexane.

This clean-up step allowed a good recovery and a low matrix effect for all kind of food examined.

Instrumental response was linear in 0.010-0.100 µg/kg range.

Specificity was evaluated on 30 different milk sample (fresh, pasteurized, UHT, bovine, ovine, liquid or powdered), 20 different infant food (milk, milk flour, cheese, dessert, yogurt) and 20 different dairy products (cheese, dessert, yogurt, ice cream): no interference was found.

Accuracy was evaluated on the same samples spiked before the extraction with AFM1 at 0.015 µg/kg for milk and infant formula and at 0.040 µg/kg for dessert and cheese. AFM1 was identified in all samples, and accuracy error was less than 5% at the interesting level.

Conclusions

The Official Method [3], based on the use of immunoaffinity chromatography columns (IAC) clean up and HPLC-FLD detection, can't to be used for infant food analysis, due to its higher detection limit. The developed method is able to detect AFM1 in infant food at concentration levels lower than 0.5 LM. Furthermore, in this case no IAC clean up is required, with a significant gain in time and in reagent costs, even for milk and cheese.

Novel Aspect

A new, faster and cheaper approach for official controls on food safety, based on a fast and very sensitive LC-MS method, AFM1 residues detection was developed and validated.

References

1. Regulation 2006/1881 EC
2. Biancardi A, Piro R, Galaverna G, Dall'Asta C., A simple and reliable liquid chromatography-tandem mass spectrometry method for the determination of aflatoxin M1 in milk, Food Additives & Contaminants: Part A, Volume 30(2013)
3. ISO 14501 2007

WP-211 / ANALYTICAL METHOD BASED ON UHPLC-ESI-MS/MS FOR QUANTIFICATION OF PLANT TOXINS IN HERBAL FOOD SUPPLEMENTS

Svetlana Malysheva - Jean-Francois Picron - Bart Huybrechts

Sciensano, Organic contaminants and additives, Brussels

Keywords: plant toxins; pyrrolizidine alkaloids; UHPLC-MS/MS; food supplements; method

Introduction

Pyrrolizidine alkaloids (PAs) and their N-oxides (PANOs), secondary metabolites of plants, are hepatotoxic to humans and animals and their ingestion in some cases might lead to death[1]. One of the sources of PA-intoxication is herbal remedies and food supplements that (un)intentionally contain material from PA producing plants. The objective of this study is to develop a simple and sensitive method for analysis of 35 PAs and PANOs in different food supplements.

Methods

The analytical method is based on ultra-high performance liquid chromatography (UHPLC) coupled to mass spectrometry (MS), as this technique can detect both tertiary PAs and their N-oxides. LC separation of the analytes is performed on a 1.7 µm BEH C18 column using alkaline mobile phase. Detection is done by a triple quadrupole MS (Xevo TQ-S, Waters) operated in positive electrospray ionization mode. The data are acquired in MRM mode.

Results

The method targets different PA types, namely senecionine-, lycopsamine, heliotrine- and monocrotaline-type compounds. The herbal dietary supplements are divided into several groups. In order to represent all the groups model matrices, namely seed powder of milk thistle, propolis, vegetable oil and herbal infusion, are selected for method development. These samples are free from PAs and PANOs or contain trace amounts. The optimization of sample preparation includes several techniques, such as a dilute-and-shoot

approach, a QuEChERS-based protocol and different sorbent materials for solid-phase extraction (SPE). The QuEChERS-method and SPE with ion-exchange mode as well as SPE with reversed-phase mode show promising results with regards to recovery; however, matrix effect is unavoidable in all cases. The ultimate aim is to select a procedure that has a minimal number of involved manipulations, assures good recovery of the target analytes and is applicable for the analysis of a variety of food supplements. The developed method will be in-house validated and limits of quantification, linearity, recovery, repeatability and reproducibility will be determined.

Conclusions

This analytical method allows for the simultaneous determination of 35 PAs and PANOs in food supplements at low levels, and can be used as a tool for collection of PA occurrence data, which is highly demanded by the European Food Safety Authority in order to perform more accurate exposure assessment.

Novel Aspect

The outlined method is able to provide information about the quality of herbal food supplements.

References

1. EFSA Panel on Contaminants in the Food Chain (CONTAM), EFSA J., 9: 2406 (2011).

For information please contact: svetlana.malysheva@sciensano.be

WP-212 / COMPARISON BETWEEN HPLC-ESI-QQQ-MS/MS AND SPME-GC-MS TECHNIQUE FOR ACRYLAMIDE DETERMINATION IN POTATOES CHIPS

Silvia Marzocchi - Sara Marziali - Maria Fiorenza Caboni

University of Bologna, Department of Agro-Food Science, Cesena

Keywords: HPLC-ESI-QTOF-MS, SPME-GC-MS, acrylamide, potatoes

Introduction

In April 2002, researchers at Swedish National Food Administration and Stockholm University published data about acrylamide content in fried and baked starch-based food [1,2]. Acrylamide formation in food is attributed to the Maillard browning reaction; because of the interaction between amino group of asparagine free amino acid and carbonyl group of reducing sugar [3].

Methods

Acrylamide in two different potatoes varieties, after frying, was determined by HPLC and GC coupled with MS. In the first case a HPLC-ESI-QTOF-MS was used, SPE extraction and chromatographic conditions were according to Calbiani et al [4]; using acrylamide- d_3 as internal standard. The method used for SPME-GC-MS was according to Ghiasvand & Hajipour [5] using a calibration curve for quantification.

Results

Fragmentation of acrylamide and acrylamide d_3 involves in NH_3 release, leading to a product ion at m/z 55 and 58 [4]. Regarding HPLC, positive electrospray proved to be the most sensitive mode for acrylamide determination [6]. Acrylamide identified in the samples was quantified using a calibration curve based on response ration of acrylamide/acrylamide d_3 . Considering GC, the quantification was made using a calibration curve, in the linear range 2-500 mg mL⁻¹ with a correlation coefficient (R^2) of 0.997. The two different analytical techniques allowed to determine different values for the same potato variety, due to the different sensitivity of the two instruments, but only with HPLC-MS were obtained results according to

literature [4,7]. In both cases the two different potatoes varieties have shown the same content difference, about 30%, probably due to the different asparagine content.

Conclusions

Both liquid and gas chromatography are able to determine acrylamide in food. This study confirmed HPLC-ESI-QTOF-MS as the most suitable chromatographic technique to determine acrylamide in potatoes sample, according to literature. At the same time, also SPME-GC-MS is able to determine acrylamide in food, but without the same sensitivity as the previous technique.

Novel Aspect

SPME-GC-MS technique was tested for acrylamide determination and quantification to find an alternative to the well-established HPLC-ESI-QTOF-MS method.

References

1. Lee M. R., Chang L. Y., Dou J, *Analytica Chimica Acta*, 582, 19-23 (2007)
2. Lagalante A. F., Felter M. A., *Journal of Agricultural and Food Chemistry*, 52, 3744-3748 (2004).
3. Lecart B., Jacquet N., Anseeuw L., Renier M., Njeumen P., Bodson B., Vanderschuren H., Richel A., *Food Chemistry*, 251, 64-68, (2018)
4. Calbiani F., Careri M., Elviri L., Mangia A., Zagnoni I., *Food Composition and Additive*, 87, 107-115 (2004).
5. Ghiasvand A. R., Hajipour S., *Talanta*, 146, 417-422 (2016).
6. Rosén J., Hellenäs K. E., *Analyst*, 127, 880-882 (2002).
7. Becalski A., Lau B. P. Y., Lewis D., Seaman S. W., *Journal of Agricultural and Food Chemistry*, 51, 802-808 (2003).

WP-213 / A MULTI-RESIDUE CONFIRMATORY METHOD FOR THE DETERMINATION OF SEDATIVES BY HPLC-MS/MS

Pavel Metalnikov - Ilya Batov - Renat Selimov - Alexander Komarov - Vyacheslav Nikulin

VGNKI, Food Safety, Moscow

Keywords: mass spectrometry, HPLC, sedatives, validation

Introduction

Sedatives are injected right before transport of the animals to reduce the stress of animals generated during handling and transport to the slaughterhouse. Thus, animal tissues have become major targets for residue analysis of sedatives and adrenergic blockers. Therefore, simple and reliable analytical methods are required to monitor these drug residues in food of animal origin. LC-MS/MS was the technique of choice due to its high specificity and sensitivity.

Methods

Samples were homogenized and fortified by addition of Internal Standard Mixture. Liquid-Liquid extraction was carried out by using acetonitrile. SPE was performed on silica gel cartridges. The HPLC-MS system consisted of HP 1290 Infinity (Agilent Technologies, USA) and 6500 QTRAP (Sciex, Canada) hybrid linear ion trap/triple quadrupole mass spectrometer. MS/MS parameters were optimized for each compound. The analysis was carried out with gradient elution.

Results

Several types of liquid-liquid extraction and SPE cartridges were tested during the study. The method was validated for muscle, liver and milk matrices. The validation study consisted of 4 individual experiments with two factors: operator and storage time before LC-MS/MS analysis. The validation levels for muscle

matrix were 1, 10, 20, 100, 200 ng/ml for azaperol, azaperone, acepromazine, haloperidol, detomidine, carazolol, xylazine, medetomidine, meperidine, triflurpromazine and fluphenazine. 5, 25, 50, 100, 250, 500 ng/ml for diazepam, metoprolol, promazine, propionylpromazine, romifidine, and chlorpromazine. Quantification was performed using individual matrix calibration curve for each series with the same concentration levels for each analyte. In order to prove the specificity and the lack of susceptibility to matrix interferences, several blank samples fortified with internal standard were additionally analysed in each series. The recovery was in the range of 45 – 74 % for all analytes in the validated concentration range.

Conclusions

Method was developed and successfully implemented for analysis of sedatives in samples of animal origin in the frames of the National Program of Monitoring of prohibited and regulated substances in food in Russian Federation.

Novel Aspect

A rapid and sensitive method for the determination of sedatives in meat, plasma and liver by HPLC-MS/MS was developed.

WP-214 / DETERMINATION OF FENPYROXIMATE FROM FOUR TYPES OF HONEY BY LC/MS/MS

Seung-Woon Myung - Jin-Mun Kim - Seung-Ho Lee

Kyonggi University, Chemistry, Suwon

Keywords: LC-MS/MS, Honey, Fenpyroximate, Pesticide

Introduction

Fenpyroximate is used to treat honey bee colonies in beekeeping farms. This not only deteriorates the quality of the honey, but also affects the health of humans who ingests it [1-2]. However, since there is no effective analytical method for the residual fenpyroximate in different types of honey, this study aims to establish an LC-MS/MS-based analytical technique after simple and rapid pretreatment with a minimum amount of sample.

Methods

5 mL of a saturated aqueous NaCl solution and 15 mL of purified water were added into 5 g of honey sample and then 10 mL of dichloromethane was added. After mixing and centrifuging, and the dichloromethane layer was then transferred to a new tube. 2 mL of the dichloromethane extract was loaded onto a SPE cartridge, and then eluted with 6 mL of dichloromethane/n-hexane (65:35%, v/v). The eluate was redissolved with 0.2 mL of acetonitrile and injected into the LC-MS system.

Results

Experiments were carried out to determine the optimal values of parameters such as sample pH, type of solid-phase cartridge, type and volume of extraction solvent, loading sample volume, type and volume of eluent for effective sample extraction, purification, and concentration. Chromatographic separation was performed using a C18 (3.0 mm i.d. × 50 mm length, 2.7 μm particle size) column under isocratic conditions of water and acetonitrile (30:70). The optimized parameters for the extraction of fenpyroximate are as follows: sample pH, 7; type of sorbent solid phase, aminopropyl; volume of extraction solvent, 20 mL; volume of loading sample, 2 mL; eluent, dichloromethane/n-hexane (65/35%, v/v); eluent volume, 6 mL. The limits of detection and quantitation for spiked honey samples were 1.56 and 5.00 μg/kg, respectively. Linear correlation coefficient (r^2) obtained from a concentration range 5.00–50.0 μg/kg was 0.9914–0.9977, while the accuracy and precision were 87.99–119.06% and 0.68–14.71 % (RSD), respectively.

Conclusions

The LOQ for the established analytical method was 5.0 µg/kg, which is 10 times lower than the EU-regulated concentration of 0.050 mg/kg. After monitoring 30 honey products, obtained from the market, using the established method, no residual fenpyroximate was detected in any sample. Compared with other studies, the novel aspects of this study are that method validation of four species of honey (acacia honey, wild honey, manuka, and chestnut honey) has been performed.

Novel Aspect

The LLE-SPE combined method does not require additional cleaning procedures, and good sensitivity can be obtained even with low-capacity cartridges.

References

1. Seenivasan S., Muraleedharan N., Int. J. Agric. Food Chem., 3, 1(2015).
2. Kiljanek T., Niewiadowska A., Semeniuk S., Gawel M., Borzecka M., Posyniak A., J. Chromatogr. A, 100, 1435(2016).

WP-215 / SUPERCRITICAL FLUID CHROMATOGRAPHY – MASS SPECTROMETRY AS A NEW TOOL FOR DETERMINATION OF SYNTHETIC FOOD DYES

Volodymyr Pauk - Martina Uhrová - Karel Lemr

RCPTM, Department of Analytical Chemistry, Palacky University in Olomouc, Olomouc

Keywords: electrospray ionization; multiple charged ions; polar analytes; supercritical phase; subcritical phase.

Introduction

Comprehensive analytical techniques are required to meet the increasingly stricter regulations on food safety. Supercritical fluid chromatography with mass spectrometric detection (SFC-MS) was assessed for analysis of food dyes. Fourteen synthetic dyes representing azo, triphenylmethane, xanthone, indigoid, and quinoline classes and covering a broad range of polarities (logD from -8 to 3.5) were selected from the list of permitted additives [1].

Methods

An Acquity UPC2SFC system coupled to a Xevo TQ-S triple quadrupole mass spectrometer with electrospray ionization (all from Waters) was used for analysis. A sharp gradient (5 to 100%) of methanolic ammonia modifier was employed on a 2-ethylpyridine column. Pseudo-SPE extraction on nylon syringe filters was utilized for real samples.

Results

Owing to the presence of multiple sulfonic groups as Na⁺ salts, multiple charged ions were observed and development of SRM transitions required careful optimization. Addition of 10 mM ammonia and 5% water into a modifier allowed separation of a wide range of analytes in a single run without using ion-pairing agents. Thus, MS compatibility, good sensitivity and wide linear range were maintained: LLOQ 4-36 ng/ml, except Brilliant black BN 255 ng/ml; ULOQ 2000 ng/ml. Mobile phase properties changed continuously from supercritical to subcritical and liquid state and did not affect the method performance, which is confirmed by validation parameters: retention time repeatability (RSD<0.6%) and peak area precision (RSD<8.9%). Brilliant blue dye was determined in a sports drink (1.9 µg/ml) and in candies (10.5 µg/g), as well as indigocarmine in candies (5.9 µg/g), which fits well within the legislative limit. The total analysis time including equilibration was 16 minutes.

Conclusions

A new SFC-MS method was developed for quantitative analysis of permitted synthetic food dyes and its essential validation parameters were evaluated. The results indicate that SFC-MS might serve as a reliable tool for analysis of both moderately polar and highly polar dyes in the context of food safety.

Novel Aspect

SFC-MS was used for analysis of ionic sulfonated dyes. A sharp modifier gradient suggested operation far from supercritical region but results were reproducible and quantitative.

Acknowledgment

The authors gratefully acknowledge the financial support from the Ministry of Education, Youth and Sports of the Czech Republic (LO1305) and Palacký University in Olomouc (IGA_PrF_2018_027).

References

1. Regulation (EC) No 1333/2008 of the European Parliament and of the Council, Official Journal of the European Union, L 354, 16-33 (2008).

WP-216 / SIMULTANEOUS SCREENING OF 6 DIFFERENT ANTIBIOTIC FAMILIES IN MEAT USING LDTD-MS/MS QUANTITATION AT 9 SECONDS PER SAMPLE

Pierre Picard - Jean Lacoursiere - Serge Auger

Phytronix Technologies inc, R&D, Quebec

Keywords: Antibiotic, Meat, screening, LDTD-MS/MS

Introduction

Food quality is a fundamental aspect of public health and is determined by innocuity associated to it. The absence of contaminants can only be validated by adequate measurement. Antibiotics in meat represent a challenge in terms of sample complexity. Until now, only a statistical sample of the total animals is tested due to the throughput. We present a new approach to screen 6 different families of antibiotics simultaneously with LDTD-MS/MS at 9 seconds per sample.

Methods

1 g of pork meat is weighed and placed in a 24 wells plate. Sample preparation is accomplished at a low temperature and reduced light. Extraction solution of MeCN:water (80:20) is added using a CERTUS FLEX to create a concentration of 1 g/mL. The samples are blended in parallel for 3 minutes and centrifuged for 3 minutes. 30 μ L of supernatant is mixed with 150 μ L of Methanol:water (75:25) and vortexed for 30 seconds. 4 μ L are deposited and dried into the analyzing Lazwell plate.

Results

Strategies have been applied to find adequate markers for thermally unstable molecules. Beta-Lactam Amoxicillin and Ampicillin are detected at the methanolic adduct formation (1). Ceftazidime undergoes a thermal fragmentation which is reproducible and the fragment containing the beta-lactam ring allows a sensitive measurement for screening purposes. The other molecules tested are detected at their pseudo molecular mass $M+H^+$. The extraction efficiency using a simple approach for the 6 types of antibiotics (Beta Lactam, Sulfonamide, Phenicol, Lincosamide, Tetracycline and Macrolides) ranges from 67 to 93%. QC samples at 50, 100 and 150% of cut-off values for each antibiotic do not overlap and exhibit sufficient signal levels for detection. Calibration curve linearity (R^2) reaches values greater than 0.993. Accuracies between 86.5% and 112.3% are observed and reproducibility ranges from 7.2 to 13.7%. Automation is optimized to achieve 96 samples in less time (12 minutes) than LDTD-MS/MS analyses (14.4 minutes) at 9 seconds per samples.

Conclusions

A complete workflow solution has been developed to extract simultaneously 6 types of antibiotic in meat. Mass spectrometer allows simultaneous specific quantitation of those compound at 9 seconds per sample using LDTD-MS/MS. Thermally unstable compounds are monitored on reproducible fragment. Sample preparation automation match the analyzing speed to enable real high throughput analysis.

Novel Aspect

Simultaneous LDTD-MS/MS screening of 6 antibiotic families in meat

References

1. Grujic S., Vasiljevic T., Lausevic M., Ast T., Rapid Communications in Mass Spectrometry, 67-74 (2008)

WP-217 / NDL-PCBS, PBDES, HBCDS LEVELS AND TEMPORAL TRENDS IN MUSSELS (MYTILUS GALLOPROVINCIALIS) FROM MIDDLE ADRIATIC SEA

Arianna Piersanti⁽¹⁾ - ***Tamara Tavoloni***⁽¹⁾ - ***Tommaso Stecconi***⁽¹⁾ - ***Arianna Stramenga***⁽¹⁾ - ***Roberta Galarini***⁽²⁾

Istituto Zooprofilattico Sperimentale dell'umbria e delle Marche "Togo Rosati", Sezione di Ancona, Ancona (1) - Istituto Zooprofilattico Sperimentale dell'umbria e delle Marche "Togo Rosati", Sede di Perugia, Perugia (2)

Keywords: PBDEs, HBCDs, NDL-PCB, mussels, temporal trends

Introduction

Mussels harvesting is a relevant commercial activity along the Adriatic Sea. Mussels are often used as sentinels for pollution monitoring, too. The aim of this study was to compare “old” contaminants like PCBs and contaminants not yet included in monitoring programs (PBDEs and HBCDs), studying temporal trends in *M. galloprovincialis* collected along the Marche Coast (2001, 2011 and 2017).

Methods

Thirteen sampling stations (9 breeding and 4 wild) along Marche Coast (Italy) were included in this study. Breeding were distributed along the Coast from North to South, while the wild ones were around Ancona. Thirty-three mussel samples were analyzed for a series of environmental pollutants: 6 PCBs, 15 PBDEs (GC-MS/MS) [1,2] and 3 HBCDs (LC-MS/MS). All data were processed using STATA software.

Results

The average concentrations of 6 PCBs were 3.75, 3.42 and 2.75 ng/g for mussels collected in 2001, 2011 and 2017, respectively and the predominant congeners were 153>138>101. No significant concentration differences were highlighted among the years. Mussels collected in the wild sites showed PCB levels significantly higher respect to the breed ones ($p < 0.001$) irrespective of the sampling period. The wild sites are closer to the coast (<1 km) and therefore subjected to higher anthropic impact. The PBDEs levels significantly decrease going from 2001, 2011 to 2017 ($p < 0.001$) [3]. The 15 PBDE medium bound decreases from 458 pg/g in 2001 to 432 and 306 pg/g, respectively, in 2011 and 2017. PBDEs contamination pattern was dominated by tetra- and penta-BDE: 47, 99 and 100 were measured in all samples and the PBDE-47 was the highest, accounting for 60% of the quantified congeners. The results were comparable with those reported in literature for European mussels [4,5]. No significant differences were observed for HBCDs among the three sampling periods. The 3 HBCDs were 273, 448 and 227 pg/g in 2001, 2011 and 2017, respectively. 2-HBCD was measured in all the mussel samples, accounting for roughly 90% of the HBCD contamination. In 20 out of 33 mussel samples, also 2-HBCD was detected [4].

Conclusions

NDL-PCB, BDE-47 and 2-HBCD were detected in 100% of the analyzed samples. PCBs did not reduce

significantly in mussels with time, but the wild sites seems more contaminated than the breeding ones. The results obtained show a PBDE decrease in the years. The BDE-47 levels were generally lower than 2-HBCD's.

Novel Aspect

This study reports for the first time the concentrations of HBCDs in Adriatic mussels, confirming the predominance of the 2-HBCD isomer.

References

1. Tavoloni T. et al, 6th International Symposium RAFA, 2013
2. Piersanti A. et al, Food Anal. Methods, 11:355–366, 2018
3. Chen D. et al, Environmental Science and Technology 45:8254-8261, 2011
4. Munsch C. et al, Chemosphere 118:329-335, 2015
5. Piersanti A. et al, Marine Pollution Bulletin 101:417–421, 2015

WP-218 / AUTOMATED SOLID PHASE MICROEXTRACTION (SPME) WITH UHPLC-MS/MS TO DETERMINE VETERINARY DRUGS IN CHICKEN MEAT

Jingcun Wu⁽¹⁾ - Abir Khaled⁽²⁾ - Feng Qin⁽¹⁾ - Janusz Pawliszyn⁽²⁾

PerkinElmer, 501 Rowntree Dairy Road, Woodbridge (1) - University of Waterloo, Waterloo (2)

Keywords: Multiclass multi-residue method, automated solid phase microextraction (SPME), UHPLC-MS/MS, veterinary drugs, chicken tissues

Introduction

Veterinary drugs are used in animal production to treat diseases and protect growth of animals. However, improper use of drugs can lead to residue violations in food products and possible health risks. Regulatory agencies around the world have established maximum residue levels or tolerance levels of veterinary drugs in foods. A common goal in drug analysis is to get acceptable results for many analytes by a cost-effective method in a single run. However, it is challenging to develop such a method for animal tissues due to the complexity of sample matrices and diversity of analytes from various classes of chemical properties. Sample preparation has always been the major bottleneck in any analytical procedure for the determination of chemical residues in food products. For veterinary drugs analysis, the most often used method has been the modified QuEChERS method with d-SPE clean-up [1-2]. However, some analytes could be lost during sample preparation due to adsorption to the sorbents. To overcome this problem, the method of direct injection of the sample extract without d-SPE cleanup has also been applied for veterinary drugs analysis recently [3]. Solid phase microextraction (SPME) is a solvent-less, effective and versatile sample preparation technique, which has been widely used in environmental and other applications including foods [4-5]. Recently, with the development of new extraction phases and devices, SPME has found increased applications in food analysis [6]. In this study, a multiclass, multi-residue method for analysis of 73 veterinary drugs (covering 13 different classes with various physical and chemical properties) in chicken muscle was developed by coupling automated SPME sample preparation with UHPLC-MS/MS. To get good analyte recovery, reduce sample matrix effects and carryover, the parameters affecting SPME performances were evaluated and optimized, such as different SPME extraction phases, extraction and desorption conditions. The method was validated according to regulatory guidelines and compared with the currently used sample preparation methods.

Methods

Chicken samples were homogenized and ground with dry ice to obtain uniform powder and kept at -20°C overnight to allow sublimation of residual dry ice and then stored at -20°C until analysis. Certified organic samples were used as matrix for method validation. 2 g of sample was weighed into a 15 mL centrifuge tube, spiked with internal standards and fortified with targeted analytes, vortex for 1 minute and leave the

sample on an agitator for 1 hour. Finally equilibrate the samples overnight in the fridge before extraction by SPME and other methods.

SPME sample preparation method: place sample tubes in the agitator for 1 hour, then add 6 mL of LCMS grade H₂O for each sample and vortex for 1 min. Transfer 1.5 mL of sample to 96 well plates. Perform static extraction for 60 min and 20 min desorption in 400 µL of desorption solvent (25:37.5:37.5 H₂O: MeCN:MeOH + 0.1% FA).

Other sample preparation methods: samples were extracted with 10 mL of extraction solution (MeCN/water, 4/1 in v/v). After centrifugation, the extract was analyzed by LC-MS/MS directly or after d-SPE clean-up with different sorbents.

UHPLC-MS/MS method: analytes were separated using an UHPLC column on a PerkinElmer UHPLC system with a gradient program with mobile phase A (water containing 0.1% formic acid) and B (MeCN containing 10% MeOH and 0.1% formic acid). Mass detection was carried out on a PerkinElmer QSight™ 220 triple quadrupole mass detector with a dual ionization source. Optimization of MS/MS parameters (choice of ions, collision energies and so on) was done by infusion of standards and use of the software. Source conditions were optimized by flow injection (FIA) method. Since the mass detector is capable of fast polarity switching, both positive and negative ionization modes were used. All instrument control, data acquisition and data processing was performed using Simplicity 3Q™ software. Calibration curves were built by standards prepared in reagents only (RO) and in sample matrix. Matrix effects were evaluated by comparison of the slopes obtained from RO and sample matrix.

Results

Mass detection method for 73 drugs was developed using time-managed-MRM algorithm of the QSight Simplicity 3Q™ software; the optimum dwell time of multiple MRM transitions was generated automatically for the targeted analytes, which saved time in method development and improved data quality and analytical performance. All calibration curves built from both the reagents only (RO) and chicken sample matrix showed good linearity (0.1X to 2X, X are the tolerance levels) with correlation coefficient (R^2) larger than 0.99 for most the analytes studied. In this study, sample matrix effects were evaluated by comparing the slopes of calibration curves obtained from RO with slopes obtained from chicken sample matrix. Sample matrix effect (%) was calculated by the percentage difference between the slopes. When the percentage difference is positive, there is a signal enhancement effect, whereas a negative value indicates signal suppression effect. It was found that sample matrix effects are significantly reduced during SPME sample preparation (only 4 out of 73 analytes have matrix effect > 20%) compared to other sample preparation method (12 out of 73 analytes have matrix effects >20%). To overcome matrix effects and reduce variations in analytical results, matrix-matched calibrations were used in all the methods for quantification of all analytes. For SPME method, the method accuracy (or relative recovery) was evaluated by spiking the analytes to the samples at three concentration levels of 0.3X, 0.9X and 2.5X, respectively and the recoveries of analytes ranged from 70% to 120 % with RSD < 20% for almost all the analytes studied, demonstrating good accuracy and precision of the method. To compare with other sample preparation methods (direct injection and d-SPE clean up), the recoveries of analytes were also evaluated by spiking the analytes to the samples at two concentration levels of 0.5X and 1X, respectively and the recoveries of analytes ranged from 70% to 120 % with RSD < 20% for most of the analytes studied when the extract was analyzed directly without clean-up. However, the recoveries for several classes of analytes are lower after d-SPE clean-up, such as for fluoroquinolones, tetracyclines and some tranquilizers. Similar results were also obtained using different sorbents such as end-capped C18 and Z-Sep+ for d-SPE in this study, which are in good agreement with the published results[1-3]. The limits of quantification (LOQs) for the drugs studied are all below the tolerance limits or the maximum residue levels (MRLs), but SPME method provided lower LOQ compared with other methods. The developed SPME-LC-MS/MS method can be applied for the fast screening and quantification of multi-class veterinary drug residues in chicken samples.

Conclusions

A multiclass, multi-residue method for analysis of 73 veterinary drugs in chicken muscle was developed by

coupling automated solid-phase microextraction (SPME) sample preparation with UHPLC-MS/MS. To get good analyte recovery, reduce sample matrix effects and carryover, the parameters affecting SPME performances were evaluated and optimized. The method was validated according to regulatory guidelines, taking into account Canadian maximum residue limits (MRLs) and US maximum tolerances for veterinary drugs in meat. A multi-day validation demonstrated that the developed method is suitable for fast and reliable quantitative analysis of 73 veterinary drugs in chicken muscle at half or below the maximum regulatory levels. Compared with other sample preparation methods such as direct injection and the frequently used modified QuEChERS method, SPME method demonstrated the following advantages: it requires less sample and solvents, easy to be automated, has less sample matrix effects and has the potential to be used for in-vivo animal sampling and fast screening by direct coupling to mass spectrometry. The selectivity and sensitivity of SPME can be tuned by using different extraction phases.

Novel Aspect

A new, fast, effective and automated SPME sample preparation method coupled with UHPLC-MS/MS for high throughput multi-class multi-residue analysis of veterinary drugs in chicken

References

1. Schneider M. J., Lehotay S. J., Lightfield A. R., *Anal. Bioanal. Chem.* 407, 4423 (2015).
2. Anumol T., Lehotay S. J., Stevens J., Zweigenbaum J., *Anal. Bioanal. Chem.* 409, 2639 (2017).
3. Lehotay S. J., Lightfield A. R., *Anal. Bioanal. Chem.* 410, 1095 (2017).
4. Souza-Silva E. A., Jiang R., Rodríguez-Lafuente A., Gionfriddo E., Pawliszyn J., *Trends Analyt. Chem.* 71, 224 (2015).
5. Xu C. H., Chen G. S., Xiong Z. H., Fan X. H., Wang X. C., Liu Y., *Trends Analyt. Chem.* 80, 12 (2016).
6. Souza-Silva E. A., Gionfriddo E., Pawliszyn J., *Trends Analyt. Chem.* 71, 236 (2015).

WP-219 / SIMULTANEOUS QUANTITATION OF FAT SOLUBLE VITAMINS IN INFANT MILK FORMULAE USING 5500 QTRAP LC-MS/MS SYSTEM

Ashley Sage⁽¹⁾ - **Jianru Stahl-Zeng**⁽²⁾ - **Alka Verma**⁽³⁾ - **Anoop Kumar**⁽³⁾ - **Manoj Pillai**⁽³⁾

SCIEX, Market Development, Warrington (1) - SCIEX, Market Development, Darmstadt (2) - SCIEX, 121 Udyog Vihar, Phase IV, Gurgaon, Haryana, India

Introduction

Fat-soluble vitamins play an important role in the stimulation of synthesis and degradation of nutrients enhancing immune function and growth performance of living organisms. In recent years there has been great emphasis on simultaneous analysis of fat-soluble vitamins. The most commonly used methods of sample preparation for the analysis of fat-soluble vitamins in foods and pharmaceuticals include sample saponification and later extraction of the vitamins from the unsaponifiable matter. Saponification process exhibits disadvantage such as oxidation, low recovery, decomposition of vitamins labile under alkaline condition. Also there are few methods reported for simultaneous determination of fat-soluble vitamins in biological fluids and food matrices, the reported methods have only enabled determination of few vitamins, which usually have been vitamins A and E and hydroxy metabolites of D2 and D3 vitamins because other vitamins are present at extremely low concentrations. The present work aims at development of highly selective and sensitive LC-MS/MS method in order to determine complete profiles of fat-soluble bioactive compounds and to design a sample treatment method without hydrolysis for the simultaneous determination of fat soluble vitamins in infant milk formulae.

Methods

Standard Stock Solution of each vitamin (A, A acetate, A palmitate, D2, D3, E, E acetate, E Succinate, K1, K2) was prepared by mass in ethanol containing 0.025% (m/v) of BHT to provide a concentration of 1 mg/mL

for all the fat soluble vitamins. Chromatographic performance was optimized using different combination of mobile phase applied to variety of C18 column to achieve optimum chromatographic resolution. Best chromatographic performance was achieved by using Phenomenex Kinetex 2.6 μm C18 (100 \times 4.6) mm column with mobile phase A: 2 mM ammonium tetra fluoro acetate and B: 100 % Methanol. . Sciex QTRAP[®] 5500+ LCMS/MS mass spectrometric system was used in multiple reaction monitoring (MRM) mode. Analysis was carried out in positive polarity with electron spray ionization mode. Multiquant[™] 3.0.2 software was used for quantitation analysis.

Preliminary Data

Protein precipitation followed by liquid-liquid extraction was used as extraction technique. Ethanol containing 0.025% of butylated hydroxytoluene was used as precipitating agent and dichloromethane was used as extraction solvent. Calibration curve with the range of 5 to 500 ppb was plotted with linear regression coefficient (R²) value greater than 0.99 for each analyte. Percent recovery ranged from 50 to 85 % for all the vitamins. Matrix effect was also calculated in terms of matrix factor and percent CV was found to be less than 15 % for each analyte. Multiquant[™] 3.0.2 analysis revealed 6 orders of linear dynamic range of quantitation. Limit of quantitation (LOQ) was achieved in nanogram level for all of the vitamins. Ion ratio of quantifier to qualifier was used for quantification.

Novel Aspect

Highly sensitive and selective method with good chromatographic resolution for simultaneous quantification of fat soluble vitamins in the concentration range of 5 ppb to 500 ppb in infant formulae without hydrolysis(saponification) has been developed.

WP-220 / MULTICLASS METHOD FOR THE QUANTITATIVE ANALYSIS OF 108 VETERINARY DRUGS IN MEAT AND MILK BY LIQUID-CHROMATOGRAPHY QUADRUPOLE-ORBITRAP MASS SPECTROMETRY

Giorgio Saluti - Irene Diamanti - Rosanna Rossi - Carolina Barola - Danilo Giusepponi - Simone Moretti - Roberta Galarini

Istituto Zooprofilattico Sperimentale dell'Umbria e delle Marche "Togo Rosati", Centre for Development and Validation of Methods, Perugia

Keywords: veterinary drugs; meat; milk; LC-Q-Orbitrap

Introduction

Multiclass approaches are more and more appreciated since they improve the cost effectiveness of monitoring programs to detect residues in food. Starting from our previous works [1; 2] in which we developed procedures for the determination of 62 antibiotics in meat and milk, the aim of this study was the widening of the method scope including 46 new veterinary drugs.

Methods

Sample preparation and instrumental conditions were those reported by Moretti et al. [1,2] with some modifications. The analysis was performed using liquid chromatography/quadrupole-orbitrap mass spectrometry. The acquisition was carried out in positive ionization mode (ESI+), except for closantel and nitroxinil (ESI-).

Results

Some slight modifications of sample treatment and chromatographic conditions were required to obtain satisfactory quantitative performances (accuracy) for the new 46 analytes, most of which were antiparasitic drugs such as avermectins and benzimidazoles. Only an example: prior to LC injection the dried residue had to be dissolved in a mixture containing ammonium acetate/acetonitrile (50:50, v/v), instead of ammonium

acetate as previously optimized [1;2]. In fact, the introduction of lipophilic drugs such as ivermectin or closantel required a certain percentage of organic solvent to avoid insufficient solubilization of these analytes and, therefore, scarce accuracy. At the end of the development, a full validation study was carried out. The observed recoveries were generally higher than 70% and precision (inter-lab reproducibility) lower than 20 % both for meat and for milk.

Conclusions

The widening of the number of analytes (from 62 to 108) included in our previously developed multiclass methods for antibiotics in food has been successfully achieved. The validation study demonstrated that the new procedure is fit for purpose [3] and therefore applicable in the official monitoring programs.

Novel Aspect

The flexibility of multiclass approaches must be experimentally verified when quantitative confirmatory methods are developed, as here demonstrated.

References

1. Moretti S., Dusi G., Giusepponi D., Pellicciotti S., Rossi R., Saluti G., Cruciani G., Galarini R. ; Journal Chromatography A, 1429, 175-188 (2016).
2. Moretti S., Cruciani G., Romanelli S., Rossi R., Saluti G., Galarini R., Journal of Mass Spectrometry, 51: 792-804 (2016).
3. Commission Decision 2002/657/EC; Official Journal of European Union, L 221, 8-36 (2002).

Acknowledgements

The authors gratefully acknowledge financial support from the Italian Health Ministry ("Determination of veterinary drugs in food: multiclass methods 2.0" IZSUM RC 022015).

WP-221 / FUNGICIDES IN RED WINES PRODUCED IN SOUTH AMERICA

Thais Vargas⁽¹⁾ - *Nathacha Salustriano*⁽¹⁾ - *Bruna Klein*⁽²⁾ - *Samantha da Silva*⁽³⁾ - *Wanderson Romão*⁽⁴⁾ - *Roger Wagner*⁽²⁾ - *Rodrigo Scherer*⁽¹⁾

University of Vila Velha, Vila Velha (1) - Federal University of Santa Maria, Santa Maria (2) - Federal University of Espírito Santo, Vitória (3) - Instituto Federal do Espírito Santo, Instituto Federal do Espírito Santo, Vila Velha (4)

Keywords: red wine; pesticides; residues; agrochemicals

Introduction

Attack of the grapevines by fungi is one of the main challenges facing producers as it compromises the production volume and quality of grapes destined for wine production. Due to these diseases, fungicides are routinely used. Thus, the objective of this work was to evaluate the presence of fungicides in 190 wines produced in southern South America to promote discussion regarding the importance of monitoring pesticide levels in wines.

Methods

A total of 190 red wines were analysed; 41 of which were Brazilian, 75 were Chilean, 30 were Uruguayan, and 44 were Argentinian. The fungicides were analysed by LC-MS/MS and GC/MS using a method validated according to the SANTE recommendations [1]. This work proposed an advancement in the quantification of dithiocarbamate in which the levels are determined in mancozeb equivalent.

Results

Of the 190 wines analysed, 76 contained fungicide residues (40%), and of those, 31 (16.3% of the total

analysed) presented values above those allowed for the grape, which would lead to them being classified as irregular. Among the 44 pesticides evaluated, dithiocarbamate was the most prevalent as it was found in 16.3% of the samples (n = 31); the next most common was azoxystrobin, which was found in 14% of the samples (n = 26), and of those, 11 were irregular. Carbendazim was found in 13% of the samples (n = 25), of which 12 were outside the limits set for grapes. Almost 20% of the samples (n = 32) contained residues of 2 or more pesticides. One of the wine samples (a blend from 2014) had five different pesticides, carbendazim, metalaxyl-M, azoxystrobin, dithiocarbamate and dimethomorph. In addition, 8 samples showed cyprodinil residues, a fungicide prohibited in grapes.

Conclusions

The main fungicides found were dithiocarbamates, azoxystrobin, carbendazim, difenoconazole, dimethomorph and cyprodinil. Of the 190 wines analysed, 40% contained at least one pesticide, and 16.3% had levels above the MRL (irregular). The present work emphasizes the importance of the creation of laws regulating pesticides in wines, which will facilitate foreign trade.

Novel Aspect

To date, this is the largest study on pesticide residues in South American wines.

References

1. European Commission, Guidance document on analytical quality control and method validation procedures for pesticides residues analysis in food and feed (2015) SANTE/11945/2015.

WP-222 / A ROBUST AND SENSITIVE METHOD FOR THE DIRECT ANALYSIS OF POLAR PESTICIDES IN FOOD AND ENVIRONMENTAL SAMPLES

Detlev Schleuder⁽¹⁾ - Ugo Chiuminatto⁽²⁾ - Wim Broer⁽³⁾

AB SCIEX Germany GmbH, Support, Darmstadt (1) - SCIEX, Support, Milano (2) - Nofalab, Laboratories, Schiedam (3)

Introduction

The prevalence of multi-residue LC-MS/MS analyses for the quantification of pesticides in food and environmental samples has been steadily increasing for many years, and they are now considered to be a minimum requirement of most laboratories working in these fields. Modern tandem quadrupoles are capable of detecting such regulated compounds at very low levels with minimal sample preparation, such as QuEChERS, thereby enabling labs to process large numbers of samples for many analytes with a fast turnaround. However, some very polar analytes which are not amenable to the extraction procedure, chromatographic method or are poor ionisers require additional single-residue methods which involve time-consuming preparation and separation and often involve derivatisation to improve detection.

Recent increase in public concern regarding the presence of glyphosate has significantly increased the requirement to analyse it and its metabolites in food, feed and the environment, so has accelerated the need for a more efficient and robust analytical method. The extraction and chromatography of these compounds is well described in the EURL-QUPPE method, but the separation is not robust in practice, so system and method maintenance are intensive. Several different HPLC or HILIC based methods have failed to address the issues of reproducibility and sensitivity, so FMOc derivatisation prior to analysis is often still employed for glyphosate, AMPA and glufosinate. Although possible to automate, this procedure is still time consuming or expensive, and is not applicable to the other polar pesticides of interest.

Methods

Nofalab is an independent sampling and testing laboratory based near Rotterdam, NL specialising in the fields of food, feed and environmental safety. The increasing pressure to provide fast, quantitative analysis

has driven them to invest in additional LC-MS/MS instrumentation and develop a new method which covers as many of these polar pesticides in a single analysis as possible. Ion chromatography has been shown to be beneficial for separation, but the need for a suppressor is detrimental to MS analysis and the problem of changing inlet systems on a heavily used mass spectrometer makes it impractical in a busy lab performing primarily reverse-phase LC. So, the final method, presented here, makes use of an IC column in a method-switching RP system with MS amenable mobile phases at pH 9. Such conditions configure glyphosate ideally for MS detection with good retention and separation of the other analytes and matrix interferences. The method meets the Sante requirements of reproducibility (<20%) and recovery (80-110%), and the LOD of the method is below 0.01 mg/kg. Excellent long-term stability and robustness were achieved throughout the validation of this method for food samples extracted by the QUPPE procedure.

Results

Where environmental samples require testing, the regulatory limits are much lower and interference from matrix more problematic in traditional analyses with a short retention time, so derivatisation is often the only option. However, since glyphosate is well retained in this new method, the potential to further develop it for direct large-volume injection was investigated in collaboration with SCIEX. By modifying the gradient conditions and optimising the injection parameters, a second method specific to environmental water samples has been developed. Although the LVI is more susceptible to changes in pH (for example, due to evaporation of mobile phase) robustness has been shown to be similarly good, and allows detection of the same suite of analytes with a LOD of <0.02 ng/l.

WP-223 / NON-TARGETED FAST SCREENING OF HUNDREDS RESIDUAL PESTICIDES ON FRUITS, EGGS AND SOIL USING AMBIENT MASS SPECTROMETRY

Jentaie Shiea

National Sun Yat-Sen University, Department of Chemistry, Kaohsiung

Keywords: non-targeted screening, pesticide, ambient mass spectrometry, miniature, high-throughput

Introduction

Thermal desorption-electrospray ionization (TD-ESI) source[1] combined with a miniature triple quadrupole mass spectrometer (Ultivo, Agilent) was applied to rapidly screen hundreds of residual pesticides on fruits and in soils and eggs without or with minimal sample pretreatment.

Methods

The TD-ESI/MS system consists of sampling probe, heated oven, electrospray device, and mass analyzer. For sampling analytes on the fruit, a stainless steel probe gently scrapes through the fruit surface for 2 cm. Since soil contains complicated components, ca. 2 μ L of the organic solvent extracts was collected on the sampling probe after rapid extraction. The probe was then inserted into a heated oven to thermally desorb the analytes on it, the analytes were subsequently carried by a nitrogen gas stream to join in an ESI plume, the analyte ions were formed by ion-molecule reactions between the charged species (H_3O^+ , CH_3OH_2^+ , $(\text{H}_2\text{O})_2^+$, etc) in the ESI plume and analytes. The analyte ions were detected by the miniature triple quadrupole mass spectrometer attached to the TD-ESI source.

Results

Without redundant extraction and separation processes, the analysis of each sample is completed in 15 s. Pesticide database was well-established by using TD-ESI coupled with a miniature triple quadrupole mass spectrometer. In addition, approximate 200 pesticides can be simultaneously detected in one analysis. The detection limit of pesticides in fruits were in sub ppm level and below the legal limit of their respective concentrations.

Conclusions

Thermal desorption-electrospray ionization/mass spectrometry is a sensitive, reproducible, qualitative, and rapidly analytical technique; it is also suitable for large-scale screening and high throughput analysis of fruits and soil to assure consumer's health.

Novel Aspect

Non-targeted screening of hundreds residual pesticides can be rapidly done by using thermal desorption-electrospray ionization/mass spectrometry.

References

1. Huang M. Z., Zhou C. C., Liu D. L., Jhang S. S, Cheng S. C., Shiea J., Anal. Chem., 85, 8956-8963 (2013)

WP-224 / THE TIME TREND OF PBDE IN BREAST MILK FROM BEIJING WITHIN THE FIRST SIX MONTHS POSTPARTUM

Shuling Song

National research center for geoanalysis, Department of organic analysis, Beijing

Keywords: polybrominated diphenyl ethers; variation trend; lactation; breast milk; Beijing

Introduction

Polybrominated diphenyl ethers (PBDEs) are flame retardants [1-2]. Numerous studies have proved that PBDEs showed the hepatotoxicity, embryotoxicity, thyroidogenic and estrogenic effects on animals [3-5]. So the residue levels of PBDEs in human are of great concern because of their potential toxicity. In this study, 7 PBDEs in breast milk from Beijing were determined, and the variation trend from 5 donors during the first six months lactation was dominated.

Methods

55 breast milks were collected from 30 donors living in Beijing over 5 years. Five of donors provided one sample per month and sampled for six months. The sample was treated with liquid-liquid extraction method (LLE). In general, 10 ml of sample was extracted by mixture of ethanol, diethyl ether and hexane, and then cleaned up by packing columns and detected by GC-NCI-MS. An external standard method was applied to qualitative and quantitative analysis.

Results

The results showed that PBDE28, PBDE47 and PBDE153 were detected in samples, and the detected rates were 23.3 % for PBDE28, 26.7% for PBDE47, and 60.8 for PBDE153, respectively. Except for PBDE28, the result was similar to the research result reported by Li J. G. [2]. The concentration of PBDE47 was ranged from 0.21 to 1.21 ng/g fat, with 0.78 ng/g fat of mean value (n=30). Compared to PBDE47, PBDE28 and PBDE153 were detected with higher average value, which were 5.40 ng/g fat for PBDE28, and were 2.73 ng/g fat for PBDE153, respectively (n=30).

The analytical results of PBDE28, PBDE47 and PBDE153 in breast milk, during the first six months lactation, indicated that the three compounds were not always detected, and the concentrations don't decrease gradually, but shaved.

Conclusions

The level of PBDEs was lower than the value reported by Li in 2010. This indicated that the residue level of PBDEs in local residents or exposure risks of local residents to PBDEs have decreased. It was well known that the diet exposure was the main exposure pathway. So we take the daily intake from diet into account for the waving change of PBDEs in breast milk during the first six month lactation.

Novel Aspect

The updated residue levels of PBDEs in breast milk from Beijing.

A time trend of PBDEs in breast milk was monitored during the first six months lactations.

The level of PBDE in breast milk may not, in some extent, showed its accumulated level in human body.

References

1. Ni K., Lu Y. L., Wang T. Y., et al. International Journal of Hygiene and Environmental Health, 216, 607-623(2013).
2. Li J. G., Yu H. F., Zhao Y. F. et al. Chemosphere, 73,182-186(2008).
3. Fair, P. A., Stavros, H. -C., Mollenhauer, M. A. M., et al. Journal of immunotoxicology, 9, 96-107(2012).
4. Walter K. M., Lin Y. P., Kass P. H., et al. BMC Veterinary Research, 13,120-131 (2017).
5. Kodavanti, P. R. S., Coburn, C. G., Moser, V. C., et al., Toxicological Science, 116, 297-312(2010).

WP-225 / SCREENING AND QUANTITATION FOR FOOD MATRICES USING SCANNING SWATH ACQUISITION

Jianru Stahl-Zeng⁽¹⁾ - ***David M. Cox***⁽²⁾ - ***Jeff Rivera Rivera***⁽²⁾ - ***Holly Lee***⁽²⁾ - ***Vanaja Raguvaran***⁽³⁾

SCIEX, MDM, Darmstadt, Germany (1) - SCIEX, Research, Concord, Canada (2) - SCIEX, Software, Concord, Canada (3)

Introduction

MSMS data is important in confidently identifying a compound in a screening assay. Obtaining an MSMS for every compound in a run can be a challenge.

SWATH acquisition can acquire MSMS that represent every precursor mass of interest at every time point. However, the deconvolution of MSMS requires that each fragment ion have a different chromatographic profile or elution time. Information dependent acquisition (IDA) will consistently trigger on the most abundant compounds, but can occasionally miss some candidates. Combining SWATH with IDA improves the identification of compounds in screening assays.

Methods

X500R, LC, pesticide kit. Acquire SWATH, IDA, and SWATH+ IDA data. Process with SCIEX OS.

Preliminary Data

SWATH acquisition uses precursor isolation windows that typically range from 5-25 Da wide. Even when two or more compounds fall within a SWATH precursor isolation window, and have similar elution times, deconvolution of the MSMS is usually possible using a PCA/PCVG based technique. This MSMS is used for library searching and subsequent compound identification.

With IDA, the precursor isolation window is 1 Da wide, which greatly reduces, but does not eliminate, the chance of two compounds being in the same window. However, if IDA does not trigger on a candidate, there is no MSMS data to confirm a compound identification with.

The amount of time spent on each cycle of acquisition is limited by how many points across an LC peak are required. With SWATH acquisition, this limit determines how wide the precursor isolation window for each SWATH experiment must be. Wider windows increase the chance of deconvolution failing, narrow windows take more time. With IDA, this limit determines how many candidates can be selected each cycle. Fewer candidates increases the chance that a candidate ion is missed, but more candidates take more time. In both cases, this trade-off follows a typical pareto principle (80/20- rule) where the majority of performance is achieved early (wide SWATH windows or small number of IDA candidates). Combining these two scan modes into one acquisition method enables a few wide SWATH windows and a small number of IDA candidates to perform well.

Novel Aspect

Data independent acquisition (SWATH) combined with IDA to improve screening assay performance.

WP-226 / LC-MS/MS VERSUS GC-MS/MS FOR HEXABROMOCYCLODODECANES (HBCDS) ANALYSIS IN FOOD

Arianna Stramenga - Tamara Tavoloni - Tommaso Stecconi - Arianna Piersanti

Istituto Zooprofilattico Sperimentale dell'Umbria e delle Marche "Togo Rosati", Sezione di Ancona, Ancona

Keywords: HBCDs, LC-MS/MS, GC-MS/MS, Food, BFRs

Introduction

HBCDs are used as additives in plastic polymers to reduce flammability; they can be released in the environment and accumulate in food chain. Technical mixture contain mainly diastereomers α , β and γ [1]. Both gas-chromatography (GC) and liquid chromatography (LC) are used for HBCDs analysis, but GC does not resolve the diastereomers. The aim of this study was to set up the best analytical method for HBCDs in food comparing GC-MS/MS and LC-MS/MS.

Methods

α , β and γ HBCD were analyzed in isotopic dilution after QuEChERS-like extraction, acidic Extrelut-NT3[®]/Si-SPE and GPC clean up. Instrumental analysis was performed in GC-MS/MS on a RTX-1614 (15m x 250 μ m x 0.10 μ m) column and PTV injection (10 μ L) [2]. A LC-MS/MS method was also developed, using a Kinetex XB-C18-100 \AA column (2.6 μ m, 100x2.10mm) and, as mobile phase, MeOH (phase A) and 2 mM CH₃COONH₄/MeOH (50:50) (phase B) [3]. The method was validated in muscle at 10, 20, 50 and 100 pg/g.

Results

Literature reports that GC does not enable the resolution of α , β and γ isomers: HBCDs can be quantified only as sum. Furthermore, the development of the GC method is extremely tricky because the 3 isomers inter-converts at T>160°C and decompose at T> 240°C or in dirty GC systems [4]. The performances for GC-MS/MS method were not satisfactory, because of the impossibility to separate the three isomers yielding in a not completely accurate quantification. During method set up, we experienced LC-MS/MS to be the method of choice for HBCDs analysis, because of the absence of thermal degradation and the ability to resolve the α , β and γ isomers. Therefore the analytical performances were systematically studied in LC-MS/MS. Accuracies, in intra-laboratory reproducibility conditions, were satisfactory for all the isomers: apparent recoveries were between 92 and 107% and precisions between 9 and 22%. Recoveries of the surrogate labelled standards were higher than 70%. No laboratory background contamination was measured.

Conclusions

The results confirmed LC-MS/MS as the method of choice for HBCDs analysis enabling to reach good performances also at 10 pg/g, which was set as method LOQ. The procedure performs well in fish products and lean terrestrial muscles. The applicability to other food matrices should be further tested.

Novel Aspect

Few data are available in Italy about HBCDs in food, therefore this method enables an extensive monitoring at very low contamination levels.

References

1. EFSA Journal, 9, 2296, 2011
2. Sales C. et al, Anal Bioanal Chem, 408:449-459, 2016
3. Abdallah M. et al, Journal of Chromatography A, 1190:333-341, 2008
4. Covaci A. et al, Environment International, 29:735-756, 2003

Acknowledgments

The authors gratefully acknowledge financial support from the Italian Health Ministry – Projects Codes:

RC0052014 and RC0052014.

WP-227 / GC-MS/MS AND LC-MS/MS APPLIED TO THE ANALYSIS OF BROMINATED FLAME RETARDANTS (PBDES AND HBCDS) IN FISH SPECIES FROM ITALIAN MARKET

Tamara Tavoloni⁽¹⁾ - *Arianna Stramenga*⁽¹⁾ - *Giampiero Scortichini*⁽²⁾ - *Tommaso Stecconi*⁽¹⁾ - *Arianna Piersanti*⁽¹⁾

Istituto Zooprofilattico Sperimentale dell'Umbria e delle Marche "Togo Rosati", Sezione Di Ancona, Ancona (1) - Istituto Zooprofilattico Sperimentale dell'Umbria e delle Marche "Togo Rosati", Sede Di Perugia, Perugia (2)

Keywords: PBDEs, HBCDs, LC-MS/MS, GC-MS/MS, Fish

Introduction

PBDEs and HBCDs are brominated flame retardants (BFRs) used as additives in plastic polymers to reduce their flammability. They can be released in the environment and accumulate in food chain. Few data are available in Italy about their environmental levels and dietary exposure [1,2]. In this study, levels of PBDEs and HBCDs were measured in marine fish samples collected in local markets.

Methods

This study involved the analysis of 15 PBDE congeners (28, 47, 49, 66, 77, 85, 99, 100, 138, 153, 154, 183, 197, 206, 209) and 3 HBCD (α -, β -, γ -HBCD) in 18 marine fish samples collected from local markets in Central Italy. The samples were analyzed either by GC-MS/MS (PBDEs) or LC-MS/MS (HBCDs) in isotopic dilution, performing a QuEChERS-like extraction followed by a two steps clean-up [3].

Results

PBDEs and HBCDs were measured above the MDL in most of the samples. Their concentrations were generally in the range of ppt, only BDE47 reached ppb level in one sample. BDE47 is the dominant congener in all species and its concentration was >LOQ in 16 out of 18 samples; in 12 samples also BDEs 100 > 99 > 49 and 154 were measured. Total PBDE contamination was estimated as sum of the 15 congeners adding 0.5 LOQ in case of not quantified. The highest PBDE levels were measured in a spiny dogfish and mackerel (Sum: 3057 and 874 pg/g respectively), while cod and hound shark showed the lowest levels (Sum: 301 and 311 pg/g respectively). The highest HBCD isomer was always the β , which was quantified in 13 out of 18 samples [4]. Only in one mackerel also α -HBCD was identified. The total HBCDs content ranged from <LOQ (cod and hound shark) to 3160 pg/g (mackerel).

Conclusions

The average concentration (N=18) was 616 pg/g for Σ PBDEs and 104 pg/g for Σ HBCDs. The highest levels for both classes were measured in spiny dogfish (lipid 7.7%), confirming that fatter species are generally the most contaminated [1] while lean ones (cod and hound shark-lipid <1%) had lower concentrations. No data are available in spiny dogfish in literature while comparable results are reported for all the other species [5,6].

Novel Aspect

Among food, fish highly contribute to PBDEs and HBCDs human intake. Few data are available in Italy about BFRs in fish, therefore these results may help in assessing background levels.

References

1. EFSA Journal, 9, 2156, 2011
2. EFSA Journal, 9, 2296, 2011

3. Piersanti A. et al, Food Analytical Methods, 11, 355-366, 2018
4. Cruz R. et al, Environment International 77, 116-131, 2015
5. Aznar-Alemany O. et al, Food and Chemical Toxicology 104, 35-47, 2017
6. Tralalon L. et al, Food and Chemical Toxicology, 104, 48-56, 2017

WP-228 / ORGANIC (PBDES,NDL-PCBS) AND INORGANIC (PB, CD, HG, AS, NI, MN)CONTAMINANTS ASSESSMENT IN WILD BOAR (SUS SCROFA) BY GC-MS/MS AND ICP-MS

Tamara Tavoloni⁽¹⁾ - ***Arianna Stramenga***⁽¹⁾ - ***Valeria Castro***⁽²⁾ - ***Martina Ciriaci***⁽¹⁾ - ***Michela Conquista***⁽³⁾ - ***Francesco Griffoni***⁽¹⁾ - ***Paolo Palombo***⁽¹⁾ - ***Giampiero Scortichini***⁽²⁾ - ***Tommaso Steccoli***⁽¹⁾ - ***Arianna Piersanti***⁽¹⁾

Istituto Zooprofilattico Sperimentale dell'Umbria e delle Marche "Togo Rosati", Sezione di Ancona, Ancona (1) - Istituto Zooprofilattico Sperimentale dell'Umbria e delle Marche "Togo Rosati", Sede di Perugia, Perugia (2) - Istituto Zooprofilattico Sperimentale dell'Umbria e delle Marche "Togo Rosati", Sezione di Tolentino, Macerata (3)

Keywords: PBDEs, NDL-PCBs, Heavy metals, wild boar, hunting.

Introduction

Wild boar is an omnivore species often used as environmental pollution sentinel. It is hunted and consumed as food, therefore may expose population to relevant levels of contaminants. This study investigate the presence of some organic (PCBs and PBDEs) and inorganic (Cd, Hg, Pb, As, Ni, Mn) pollutants in wild boars caught in Central Italy.

Methods

Samples were caught in a rural area close to the Appennini mountain in the Macerata district (Marche). Muscle, fat and liver of 44 wild boars hunted in 2017/18 season were selected. The samples were grouped by organ type, sex and weight to yield 21 laboratory pools. NDL-PCBs and PBDEs were analyzed by GC-MS/MS [1,2], Pb, Cd, As, Ni and Mn by ICP-MS and Hg by AAS (FIMS).

Results

PCBs, Pb, Cd and Mn were above LOQ in all the selected tissues, PBDEs only in muscle, As was not detected in almost all the samples while traces of Ni and Hg were measure only in livers. PCBs were all below MRLs [3] and the highest levels measured in muscle (mean 13 ng/g fat), in liver mean PCBs was 0.5 ng/g; congener 153 was dominant, followed by 138 and 180. The concentration increases with weight in males while decreases in females, where the levels are generally lower, probably because they can transfer contaminants to their progeny [4]. Also for PBDEs, muscle was more contaminated than liver where almost all the BDEs were below LOQ. In muscle dominant congeners were 47> 99> 100. 15PBDEs in muscle (mean 848 pg/g) was generally higher than the ones reported from other authors [5]. Cd in muscle (mean 0.028 mg/Kg) was lower than MRL (0.050 mg/Kg) except for two samples (0.057 and 0.053 mg/Kg). None of the livers exceeded MRL (mean 0.25 mg/Kg) [3]. The results are comparable with previous studies [6, 7]. Pb values fell between 0.025 and 6.99 mg/Kg in muscle while in liver the mean level was 0.04 mg/Kg. Six out of 7 muscles exceeded Pb MRL; this could be related to Pb dispersion in muscle due to bullet fragmentation. Mean levels of Mn were 0.482 (muscle) and 1.774 mg/Kg (liver).

Conclusions

Muscle resulted to be the most contaminated tissue respect to liver and fat for all the considered contaminants. Pb residues in muscle were very high as a result of fragmentation of bullets used in hunting.

Novel Aspect

Few data are available about contamination of wild fauna in Central Italy, therefore these results may help

to start assessing background levels in Italy.

References

1. Tavoloni T. et al, 6th International Symposium RAFA, 2013
2. Piersanti A. et al, Food Analytical Methods, 11:355–366, 2018
3. European Commission regulation 1881/2006 and amendments
4. Neila C. et al, Ecotoxicology and Environmental Safety, 140-24:29, 2017
5. Zacs D. et al, Science of the Total Environment, 621:308-316, 2018
6. Chiari M. et al, Food Additives and Contaminants: Part A, 32:1867-1874, 2015
7. Amici A. et al, Journal of Animal Science, 11:353-362, 2012

Acknowledgments

The authors gratefully acknowledge financial support from the Italian Health Ministry - Projects Codes: RC0022016

WP-229 / SCREENING AND QUANTITATIVE ANALYSIS OF NON INTENTIONALLY ADDED SUBSTANCES ELUTED FROM POLYETHYLENE TEREPHTHALATE BOTTLE BY HIGH-RESOLUTION MASS SPECTROMETRY

Atsushi Yamamoto⁽¹⁾ - *Taro Murakami*⁽²⁾ - *Motoshiro Shizuma*⁽³⁾ - *Eri Kishi*⁽²⁾ - *Asako Ozaki*⁽²⁾

Tottori University of Environmental Studies, Faculty of Environment, Tottori (1) - Osaka Institute of Public Health, Food Chemistry 2, Osaka (2) - Osaka Research Institute Industrial Science and Technology, Biomaterial, Osaka (3)

Keywords: Food safety, NIAS, high-resolution mass spectrometry, fragmentation

Introduction

Food packaging materials can contain a wide range of known and unknown substances [1]. Chemicals such as impurities of ingredient, by-products, and degradation products are also plausible. These chemicals are referred to as non-intentionally added substances (NIAS). In this study, NIAS from polyethylene terephthalate (PET) bottles were examined by liquid chromatography high-resolution mass spectrometry (LC/HRMS).

Methods

To understand migration of NIAS to bottled beverages, non-volatile polar NIAS was assigned as a main objective in this study. An LC/HRMS system comprising ExionLC AD and X500R (SCIEX) was used for non-target analysis. X500R mass spectrometer, which had a mass accuracy of 2 ppm, was used to acquire the HR mass spectra with a mass resolution of 30,000. Material tests and long storage tests of PET bottles obtained in Japan, were carried out.

Results

Several compound series with a constant interval of m/z 192.042 were present. This value corresponded to $C_{10}H_8O_4$, a monomer unit of PET. The mass spectrum includes much information such as fragmentation, exact mass, neutral loss, adduct ion formation, and isotopic pattern. Each element sticks to fundamentals and overlapped one another, so that the mass spectrum looks complex. It is cumbersome but usually possible to interpret the mass spectrum to determine a structure. The congener with the highest abundance in the material tests was cyclic PET trimer. Monomers and oligomers whose one terminal was carboxy group were also identified by the interpretation of mass spectra and confirmed by reference standards. One-year long storage tests with simulated food solvent showed a gradual release of these substances to the contents. PET bottles obtained from 8 manufacturers showed roughly similar migration profiles.

Conclusions

Unlike other plastic materials and articles, PET bottles obtained in Japan included less additives. In the examination by non-target analysis using LC/HRMS, PET congeners were the main migrants from PET bottles. One-year long storage tests with simulated food solvent showed a gradual release of these substances to the contents. The highest concentration of 4-[(2-hydroxyethoxy)carbonyl]benzoic acid was 4.6 ng/mL and detected in one-year long storage test sample.

Novel Aspect

NIAS from PET bottles was examined by non-target analysis with LC/HRMS. Several PET congeners were identified from HRMS data and confirmed by reference standards.

References

1. Nerin C., Alfaro P., Aznar M., Domeño C. *Anal. Chim. Acta* 775, 14–24 (2013).

WP-230 / DEVELOPMENT OF A RAPID METHOD FOR ANALYSIS OF HETEROCYCLIC AROMATIC AMINES, BY QUECHERS EXTRACTION AND UHPLC-APCI-MS/MS ASSAY

Sylvie Chevolleau - Jade François - Alyssa Bouville - Laurent Debrauwer

INRA /Toxalim, Metatoul-AXIOM Platform, Toulouse, France

Keywords: Heterocyclic Aromatic Amines, QuEChERS, UHPLC-APCI-MS/MS, meat, cooking

Introduction

Heterocyclic Aromatic Amines (HAs) are neo-formed compounds during the cooking of meats and are known or suspected to be mutagenic and carcinogenic for rodents and man. More than 20 HAs have been classified according to their structures and physico-chemical properties. They have been characterized at concentrations in the ng/g (ppb) range in most meat cooked under normal time and temperature conditions, but the analysis procedure is still long and tedious.

Methods

We present the development of a methodology for the quantification of 16 HAs by UHPLC-APCI-MS/MS after QuEChERS extraction and SPE purification of beef meat samples. Protocols, using an internal standard, were optimized in order to reach quick analysis, a good repeatability and a maximum recovery yield for each HAs. Quantification was achieved by UHPLC-APCI-MS/MS (MRM mode) on a triple quadrupole instrument by monitoring two specific transitions for each HAs.

Results

The objective of this work was to develop a simple and fast method to overcome the traditional time-consuming extraction method of HAs using diatomaceous earth and large volumes of solvent. For that, the performance of the QuEChERS) method combined with UHPLC-APCI-MS/MS was investigated. Due to, the particular physicochemical properties of the HAs, as well as the fatty nature of the beef matrix, the classical QuEChERS method was modified, by carrying out a solid phase extraction(SPE) in mixed mode combining a polymer adsorbent and a cation exchange, for the purification of the extract in order to minimise matrix effects and signal extinction.

This method takes advantage of the rapidity, specificity and sensibility of both QuEChERS and UHPLC-APCI-MS/MS. Its performances will be presented in terms of selectivity, sensitivity, linearity, carry-over, recovery, matrix effect, repeatability, accuracy and intermediate precision in doped beef meat extracts (15% fat), according to the European Medicines Agency (EMA) validation guidelines.

Conclusions

An original methodology has been developed to quickly and easily assess 16 HAs in complex beef meat

matrices by UHPLC-APCI-MS/MS. After validation, this method will be applied to implement a strategy allowing a reasoned choice of curative ingredients (marinades ...) to mitigate the formation of HAs during the cooking of meat products and provide possible recommendations to minimize the formation of HAs during meat cooking and to reduce consumer exposure towards HAs.

Novel Aspect

Easy rapid and repeatable extraction and purification method suitable for complex and fatty meat matrix combined with a sensitive and specific UHPLC-APCI-MS/MS assay method.

WP-231 / IN VIVO EFFICACY OF A YCW-BASED PRODUCT IN REDUCING ORAL BIOAVAILABILITY OF ZEARELENONE IN PIGS

Vito D'Ascanio⁽¹⁾ - *Amelie Catteuw*⁽²⁾ - *Elisa Santovito*⁽³⁾ - *Virginie Marquis*⁽³⁾ - *Siska Croubels*⁽²⁾ - *Donato Greco*⁽³⁾ - *Antonio F. Logrieco*⁽¹⁾ - *Giuseppina Avantaggiato*⁽¹⁾

Institute of Sciences of Food Production (ISPA), National Research Council (CNR), Via Giovanni Amendola, 122/O, 70126 Bari, Italy (1) - Department of Pharmacology, Toxicology and Biochemistry, Faculty of Veterinary Medicine, Ghent University, Salisburylaan 133, 9820 Merelbeke, Belgium (2) - Phileo Lesaffre Animal Care, Marcq-en-Baroeul, France (3)

Keywords: Yeast cell wall; efficacy; zearalenone; toxicokinetic; piglets; biomarker.

Introduction

The use of yeast cell wall (YCW) derivatives as mycotoxin binders (MB) can be considered a promising approach to protect livestock against the harmful effects of mycotoxin exposure. Several in vitro studies confirmed the potential of YCW in binding mycotoxins. As in vitro studies do not always predict in vivo results, the European Food Safety Authority (EFSA) guidelines state in vivo trials are needed. This work was aimed to determine the in vivo efficacy of a YCW, obtained from Lesaffre International, in reducing oral bioavailability of zearalenone (ZEN) in piglets.

Methods

In order to investigate the in vivo efficacy of a YCW in reducing ZEN bioavailability, a toxicokinetic study focusing on the plasma concentration-time profiles of ZEN was performed. The study was conducted on 32 healthy 5-week-old piglets, sexes equally divided. The pigs were randomly allocated in 2 groups of 16 animals, a control group (C) (treated with a single oral bolus of 0.3 mg/kg BW ZEN) and a treatment group (T) (treated with a single oral bolus of 0.3 mg/kg BW ZEN in combination with the detoxifier).

Results

Plasma concentrations of ZEN were quantified with a validated LC-MS/MS method. Limit of detection (LOD) and limit of quantification (LOQ) were 0.1 and 0.5 ng/mL, respectively. As plasma concentrations of ZEN after oral administration were below the LOQ, ZEN-glucuronide was quantified as plasma biomarker for ZEN. Non-compartmental toxicokinetic modeling was performed using Phoenix 64 (Certara, USA). The effect of the mycotoxin detoxifier on the oral absorption of the mycotoxin was evaluated by statistical comparison of toxicokinetic parameters between the C and T groups, with special emphasis on area under the curve from time zero to infinite ($AUC_{0 \rightarrow \infty}$), maximal plasma concentration (C_{max}), time at maximal plasma concentration (T_{max}) and relative oral bioavailability (F).

The mycotoxin detoxifier significantly reduced oral ZEN absorption in piglets (p value < 0.05).

Conclusions

A YCW based product, produced at lab and industrial level, was selected after showing high efficacy in adsorbing ZEN in vitro. In vivo efficacy testing of the YCW based mycotoxin detoxifier revealed a significant

reduction of ZEN oral bioavailability in piglets.

Novel Aspect

In vivo efficacy testing of a safe YCW based product to be proposed as detoxifier for reducing ZEN exposure, and consequently estrogenic effects of ZEN, in farm animals.

Acknowledgments

The present work has received funding by the European Union's Horizon2020 Research and innovation programme under Grant Agreement No. 678781 (MycoKey).

WP-232 / EXPLOITING THE POTENTIAL OF HR-MS FOR DECIPHERING MYCOTOXIN MODIFICATION IN MICROPROPAGATED WHEAT: THE CASE OF DEOXYNIVALENOL

Tito Damiani⁽¹⁾ - *Laura Righetti*⁽¹⁾ - *Enrico Rolli*⁽²⁾ - *Michele Suman*⁽³⁾ - *Renato Bruni*⁽¹⁾ - *Chiara Dall'Asta*⁽¹⁾

University of Parma, Department of Food and Durg, Parma, Italy (1) - University of Parma, Department of Biosciences, Parma, Italy (2) - Barilla G. R. F. Ili SpA, Advanced Laboratory Research, Parma, Italy (3)

Keywords: deoxynivalenol, wheat, metabolism, HRMS, untargeted

Introduction

Deoxynivalenol (DON) is a mycotoxin produced by *Fusarium* sp. Plants mechanisms to counteract mycotoxins include phase I and II reaction. The resulting modified forms have different structures and chemical formulas, making difficult their detection by routine analytical approaches. In this study, liquid chromatography coupled to high resolution mass spectrometry (LC-HRMS) was employed to identify DON metabolites produced by in-vitro micropropagated wheat plantlets.

Methods

Two wheat varieties (Kofa and Svevo) were grown in a DON-spiked medium. After 14 days, media and plant samples underwent to LC-HRMS analysis. DON absorption and translocation to the aerial part were evaluated. Phase I and II metabolites were hypothesized, respectively, by Sites of Metabolism (SOM) prediction and based on the known plant detoxification pathways. Their identification was achieved using accurate mass, isotopic pattern and HRMS/MS spectra.

Results

At the end of the experiment, in vitro plantlets showed no visible damages. After 14 days, a reduction of 10% and 30% of DON was observed in the growing media of Kofa and Svevo, respectively. The more efficient DON absorption of Svevo is in accordance with previous results [1].

Quantification of DON was carried on different plant organs to evaluate the mycotoxin absorption and translocation. In the plant material, an extensive conversion of the DON into its conjugate forms, among them DON-3-Glucoside (DON-3-Glu) and other analogues, was observed.

Untargeted profiling was carried out to identify the DON biotransformation products. Metabolites annotation was performed as already proposed by other authors [2]. Full identification was achieved when commercial standards were available. For the remaining DON metabolites, they were putatively identified considering accurate mass, isotopic pattern distribution and the comparison between experimental and in silico predicted fragmentation pattern (<http://cfmid.wishartlab.com>).

Conclusions

Although the conditions applied do not represent the complexity of an actual environment, in vitro wheat culture has proved to be suitable for the investigation of mycotoxin metabolism in plant.

Svevo showed a larger DON absorption from the growing medium. However, no considerable differences

were observed in the DON-3-Glu conversion ability, between the two varieties.

A considerable number of the hypothesized metabolites was identified through the untargeted LC-HRMS analysis.

Novel Aspect

The integration of in-vitro plant cultures and HR-MS is an effective workflow for the elucidation of the metabolic fate of mycotoxins and other xenobiotics.

References

1. Righetti L., Rolli E., Galaverna G., Suman M., Bruni R., Dall'Asta C., PLoS ONE, 12(11) (2017).
2. Kluger B., Bueschl C., Lemmens M., Michlmayr H., Malachova A., Koutnik A., Maloku I., Berthiller F., Adam G., Krska R., Schuhmacher R., PLoS ONE, 10(3) (2015).

WP-233 / SENSITIVE MASS SPECTROMETRIC METHOD FOR THE DETECTION OF PALYTOXINS IN SHELLFISH USING LITHIUM CATIONIZATION

Arjen Gerssen - Mirjam D. Klijnstra

RIKILT, Wageningen University and Research, Wageningen, Netherlands

Keywords: Cationization, palytoxins, food safety

Introduction

Palytoxins (PITXs) are produced by marine dinoflagellates, cyanobacteria and certain corals and is a neurotoxic compound that can end-up in shellfish. PITXs are within the classes of marine biotoxins the largest non-protein molecules with m/z 2700. To protect shellfish consumers EFSA established an acute reference dose of 0.2 $\mu\text{g}/\text{kg}$ b. w. which can be translated to 30 μg PITX-eq/kg shellfish (400 g shellfish by 60 kg person). Sensitive methods for analysis are lacking.

Methods

Shellfish homogenate was extracted in triplicate with 50:50 v/v water/methanol. Subsequently extracts were purified with SPE clean-up. For the LC chromatographic separation a Phenomenex Kinetex (2.1 x 100 x 1.7) column was used. Mobile phases consisted out of 0.25 mM LiI and 0.31 mM formic acid in water and 0.25 mM LiI and 0.31 mM formic acid in 90% acetonitrile. The LC effluent was directed into an Q-Orbitrap-MS or tandem quadrupole MS both equipped with ESI operating ESI+.

Results

When analysed by mass spectrometry complex spectra are obtained with multiple adducts, charges states and in source water losses. To gain sensitivity there was attempted to simplify the spectrum by the formation of lithium adducts. Without the addition of lithium in the mobile phase the most abundant ion present is less than 30% of the sum of all ions formed. However with the addition of lithium the main adducts formed were $[\text{M}+\text{H}+2\text{Li}]^{3+}$ and $[\text{M}+\text{H}+\text{Li}]^{2+}$ where the abundance of $[\text{M}+\text{H}+2\text{Li}]^{3+}$ was around 60% of the total abundance of all ions formed. In MS/MS the most abundant fragments are mainly due to losses of H_2O . Also more structural fragments were obtained, the most sensitive transition is $[\text{M}+\text{H}+2\text{Li}]^{3+}$ m/z 897.8 > 1215.7. When measured with MS/MS the method has an LOQ of 1 ng/mL for a standard solution. Subsequently the method was successfully validated for PITX in shellfish at levels of 10, 30 and 60 μg PITX/kg.

Conclusions

A sensitive method based on cationization with lithium for the detection of PITXs in shellfish was developed. Subsequently the method was validated with good recoveries (70-120%) and repeatability RSDs

below 20%. The method described is a novel method for the detection of intact palytoxins in shellfish which can measure these compounds at the relevant levels ($\mu\text{g}/\text{kg}$ range).

Novel Aspect

Application of cationization in the field of marine biotoxin research for the development of an intact palytoxin LC-Q-Orbitrap and LC-MS/MS method.

For information please contact: arjen.gerssen@wur.nl

WP-234 / DEVELOPMENT OF A SIMPLE MULTI-RESIDUE DETERMINATION METHOD OF 80 VETERINARY DRUGS IN FISH TISSUES BY LIQUID CHROMATOGRAPHY COUPLED TO QUADRUPOLE ORBITRAP MASS SPECTROMETRY

Xin Luo - Zhixu Tang - Zengyuan Niu - Xiwen Ye - Fei Zhao

Qingdao Customs, Inspection and Quarantine Technical Center, Qingdao, China

Keywords: Veterinary drugs; Multi-residue; Fish tissues; Extraction optimization; Orbitrap

Introduction

There is no doubt that veterinary drugs are very important for aquaculture. But inappropriate and commonly subtherapeutic use of veterinary drugs may result in the presence of drug residues in the fish tissues even food products. Most of the published multi-residues detection methods commonly employed the LC-triple-quadrupole tandem MS techniques, which requires extensive compound-dependent parameter optimization. And due to the great differences in the physicochemical properties of the various targets, we need an appropriate extraction procedure to ensure that all compounds can be extracted effectively. On the other hand, protein, fat and many other interferences are abundant in fish. There is an effective clean-up procedure still needed to eliminate background matrix interferences.

Methods

A simple, rapid and sensitive multi-residue analytical method was developed and validated for 80 veterinary drugs in fish tissues using ultrahigh performance liquid chromatography-Orbitrap high resolution mass spectrometry (LC-HRMS). The analytes belong to 12 different families include benzimidazoles, β -lactams, lincosamides, macrolides, nitromidazoles, quinolones, sulfonamides and trimethoprim, tetracyclines, triphenylmethane dyes, amphenicols, nonsteroidal estrogens and steroid hormones. The sample preparation was optimized base on QuEChERS (quick, easy, cheap, effective, rugged and safe) procedure. A very simple and sufficient preparation procedure without salting-out and complex clean-up process was studied. It had been proved that water in the extract was helpful for extracting hydrophilic compounds and precipitating the lipids during the subsequent cleaning process. In addition, an appropriate percent of methanol was necessary to some analytes.

Results

Finally, a mixture of acetonitrile, methanol and water (3:1:1, v/v/v) which include 1% acetic acid and 10 mM ethylenediamine tetraacetic acid disodium salt 2-hydrate was selected as the extraction solvent, and the clean-up step consisted of a low temperature procedure and two times of high-speed centrifugation to deproteinize and remove lipids. The detection and quantification of all compounds were performed by ultrahigh performance liquid chromatography coupled with electrospray ionization quadrupole Orbitrap high resolution mass spectrometry in positive and negative ion mode. This methodology was validated according to the Commission Decision 2002/657/EC and SANTE/11945/2015. The recoveries ranged from 60.74%-109.85% with relative standard deviations (RSDs) < 20%. The limits of quantification (LOQs) were 0.25 to 25 $\mu\text{g}/\text{kg}$, for the analytes which the MRL or MRPL had been established in fish tissue, the LOQs

were all lower than their own legal tolerances. The values of decision limit ($CC\alpha$) and detection capability ($CC\beta$) were in the range of 1.91-1001.13 $\mu\text{g}/\text{kg}$ and 3.52-1002.26 $\mu\text{g}/\text{kg}$, respectively. This validated method has been successfully applied on the determination of veterinary drugs in real commercial oplegnathus punctatus samples.

Conclusions

A very simple, rapid and efficient multi-residue and multi-class analytical method for the simultaneous determination of 80 veterinary drugs in fish tissue by ultrahigh performance liquid chromatography coupled with quadrupole Orbitrap high resolution mass spectrometry has been developed. To the best of our knowledge the proposed method is the simplest way reported in literature for multi-residue and multi-class determination in fish. The method was validated with linearity, recovery, precision, sensitivity and tested the $CC\alpha$ and $CC\beta$ for each analyte. The satisfactory performances were achieved mainly depended on the high efficiency of extraction solution and the excellent selectivity and sensitivity of Orbitrap. This efficient and rapid method is a good way to screening for large batches of fish samples.

WP-235 / COMPARISON OF DIFFERENT STRATEGIES FOR THE DETERMINATION OF VETERINARY DRUG RESIDUES IN FOOD BY UPLC-MS/MS USING ORBITRAP ELITE

Pavel Metalnikov⁽¹⁾ - *Renat Selimov*⁽¹⁾ - *Alexey Efimov*⁽¹⁾ - *Ilya Batov*⁽¹⁾ - *Ayshat Dzhatdoeva*⁽¹⁾ - *Alexander Komarov*⁽²⁾

VGNKI, Food Safety, Moscow, Russian Federation (1) - VGNKI, Moscow, Russian Federation (2)

Keywords: Mass Spectrometry, Orbitrap, Veterinary Drugs

Introduction

Analysis of veterinary drug residues in food matrices usually is performed for a specific group of drugs. Now, with the advent of high resolution mass spectrometry, it's getting possible to use accurate masses, measured during analysis, for unambiguous identification of potential contaminants, and so-called "general unknowns" as well. The aim of the work was to compare different strategies for identification of analytes of different drug groups in large scale.

Methods

Samples of food matrices (milk or meat) were originally homogenized, fortified with veterinary drug mixture solutions, and subjected to liquid-liquid extraction. Then, Solid Phase Extraction (SPE) was applied. Analysis was done on UPLC-Orbitrap system (Ultimate 3000 – Orbitrap Elite (Thermo Scientific, USA)) in gradient mode.

Results

Orbitrap is a complex mass spectrometer, measuring masses of peaks with accuracy about 1 ppm and providing MS/MS data of high quality. Veterinary drugs were identified either by determination of their elemental composition based on accurate mass determination from spectra, detected in electrostatic ion trap (FT mode), or by fragmentation spectra, registered in linear ion trap (IT mode). Fragmentation spectra matched to in-house built Spectra Library or to other available data bases.

Conclusions

Despite accurate mass determination is enough for compound identification, in some cases (*e.g.* overlapping peaks) confirmation by fragmentation spectra is important.

Novel Aspect

New strategies for veterinary drug residues determination in food matrices were developed.

WP-236 / UHPLC-HRMS ANALYSIS OF THEOBROMINE IN COSTA RICAN THEOBROMA CACAO

Katarina Mladenovic - Yuriko Root - Dil Ramanathan

Kean University, Union, New Jersey, USA

Keywords: High Resolution Mass Spectrometry, Costa Rican Cacao, UHPLC-HRMS, Theobromine

Introduction

Costa Rican Cacao has been a staple within Costa Rican culture for many centuries. Within the past 10 years, Costa Rican cacao has made a return for chocolate enthusiasts, helping to provide antioxidants to humans and satisfy sweet cravings. One of the major components of chocolate is Theobromine. Theobromine is a bitter alkaloid/flavonoid beneficial in the treatment of hypertension, arteriosclerosis, and angina pectoris.

Methods

The samples were centrifuged and vortexed numerous times to give a purified product. A three-minute gradient method with a flow rate of 300 $\mu\text{L}/\text{min}$ was developed on the UHPLC-HRMS using HPLC-grade water and acetonitrile. Quantitative analysis on Theobromine was completed on the Thermo Scientific UHPLC and LTQ Orbitrap Discovery equipped with an ESI ion source. 2-Hydroxyethyl theophylline was used as the internal standard.

Results

In a previous study from our laboratory, Headspace GC-MS analysis was used to identify flavonoids from Cocoa Beans. Theobromine was the major component found. To give the precision of the theobromine extraction process, the recovery analysis was 86%. Raw, unroasted, roasted, and 100% cacao was analyzed in concentrations of ng/mL. Raw cocoa had a concentration of 586 ng/mL and unroasted cocoa had a concentration of 747 ng/mL. Roasted cocoa had a concentration of 521 ng/mL and 100% cacao had a concentration of 1018 ng/mL. Roasted cocoa had the lowest concentration of theobromine.

Unroasted cocoa had a higher concentration of theobromine prior to roasting cacao. Theobromine could have been lost based on the roasting process or choice of roasting. The 100% cacao chocolate bar had the highest amount of theobromine. Additional theobromine existed and contributed to a more bitter and rich cocoa flavor. Theobromine levels in popular commercially available chocolates such as Lindt, Hershey, Mars, etc will be discussed.

Conclusions

Theobromine levels in popular commercially available chocolates such as Lindt, Hershey, Mars, etc. were determined by UHPLC-HRMS and compared between raw, unroasted, and roasted chocolates.

Novel Aspect

A sensitive and selective UHPLC-HRMS method to quantify Theobromine levels in chocolates.

WP-237 / LC-MS/MS METHOD TO EVALUATE THE PRESENCE OF EXOGENOUS GLYCEROL IN WINES

Loris Tonidandel - Sergio Moser - Debora Trainotti - Roberto Larcher

Edmund Mach Foundation, CTT, San Michele all'Adige, Italy

Keywords: LC/MS/MS, glycerol, wine, cyclic diglycerols, 3-methoxypropane-1,2-diol

Introduction

Glycerol is a natural constituent of wine being produced by yeast during fermentation and its concentration ranges between 5 and 11 g/L depending on yeast strain. It plays an important role on the sweet taste of

wine [1]. The addition of synthetic glycerol to wines in order to increase this parameter is not allowed by the European oenological legislation [2]. Actually, the strategy to investigate the presence of exogenous glycerol in wine is based on the detection of impurities produced by the industrial processes on its production [3], which are absent in wines. Synthetic glycerol contains the follows impurities: 3-methoxypropane-1,2-diol (3-MPD) and cyclic diglycerols (CycDs) formed during the purification of the rough industrial product conducted through a distillation process. Currently, the determination of these impurities are carried out by an GC/MS method after an extraction step in ethyl ether. As an alternative to this "old" method has been proposed an innovative approach based on liquid chromatography coupled to the mass spectrometry (LC-MS/MS).

Methods

The chemical analysis of 3-methoxypropane-1,2-diol (3-MPD) and cyclic diglycerols (CycDs) has been performed using an Acquity Liquid Chromatographer (UPLC Waters Corporation, Milford, MA, USA) coupled to an Xevo TQ Mass Spectrometer (Waters Corporation, Milford; MA; USA). Chromatographic separation has been performed with an C18 HSS T3 column (2.1 mm x 100 mm, 1.8 µm); flow rate 0.45 ml/min; Eluent A, Water + NH₄Ac (5mM); Eluent B, Methanol. Mass spectrometer was equipped with an electrospray ion source operating in positive ion mode; capillary voltage 1.5 kV; nitrogen gas flow, 1000 L/h; source temperature, 150°C. Acquisition were carried out in MRM (multiple reaction monitoring).

Results

Gas chromatography mass spectrometry (GC-MS) is the analytical method adopted by the OIV to detect the fraudulent addition of glycerol by measure these two contaminants [4]. To the best of our knowledge, this is the only validated method to detect and quantify 3-MPD and CycDs in wine. However, the application of the GC-MS method shows serious difficulties: sample preparation and time analysis (runtime about 42 min); dirty injections with early consumption of the column (particularly in the case of sweet wines); mass spectrometer ion source must be cleaned very often to improve the required sensitivity. In order to overcome these limitations a simple LC-MS/MS method has been developed.

Conclusions

The proposed method could be a useful tool to rapidly screening wine suspected of fraudulent glycerol addition. Keeping GC-MS analysis to further confirm positive results.

Novel Aspect

High throughput LC-MS/MS method respect the official method (using GC/MS) proposed by OIV.

References

1. Noble A. C., Bursick G. F., Am. J. Enol. Vitic. 1984, 35, 110–112.
2. Commission Regulation (EC) No 606/2009 of 10 July 2009 laying down certain detailed rules for implementing Council Regulation (EC) No 479/2008 as regards the categories of grapevine products, oenological practices and the applicable restrictions.
3. Fahl C., Wittkowski R., Lofthouse J., Hird S., Brereton P., Versini G., Lees M., Guillou C. . J AOAC Int. 2004, 87, 1179-1188.
4. OIV-MA-AS315-15. Determination of 3-methoxypropane-1,2-diol and cyclic diglycerols (by-products of technical glycerol) in wine by GC-MS - description of the method and collaborative study - (Resolution Oeno 11/2007)

WP-238 / APPROACHES TO MINIMIZE AND CONTROL EXTERNAL CONTAMINATION DURING SAMPLE HANDLING AND ANALYSIS OF UBIQUITOUS ENVIRONMENTAL PHENOLIC COMPOUNDS IN FOOD SAMPLES

Aida Zuberovic Muratovic - Urszula Rybinska - Karl-Erik Hellenäs

National Food Agency, Department of Chemistry / Division of Science, Uppsala, Sweden

Keywords: food analysis, external sample contamination, ubiquitous environmental phenolic compounds, UPLC-ESI-MS/MS

Introduction

Phenolic compounds belong to a group of substances that might be of synthetic or natural origin. Synthetic phenolic compounds are present in our environment in pharmaceuticals, preservatives, personal-care products (e.g. synthetic fragrances and surfactants), as pesticides and miscellaneous industrial chemicals and byproducts (e.g. from plastic industry) [1,2].

As such, they can be detected in food and beverages in which they might originate from additives, food processing and packaging materials or as contaminants from the environment including the analyst. Many of the synthetic phenolic compounds are suspected to have the ability to affect the human body acting as estrogen mimics and are collectively termed endocrine disrupting compounds (EDCs). EDCs are substances that “interfere with the synthesis, secretion, transport, binding, action or elimination of natural hormones in the body that are responsible for the development, behavior, fertility and maintenance of homeostasis (normal cell metabolism)”[3]. EDCs are still not regulated at the EU level, but the work is ongoing to develop scientific criteria and properties for regulation. This causes intense discussion over these substances’ toxicity and mutagenicity which is reflected in ever growing interest in the question of their presence in foods. Thus, analytical methods for reliable quantitative determination of ubiquitous environmental phenolic compounds in food are desirable, but the analysis represents a substantial challenge due to the overwhelming risk of unintentional sample contamination from the laboratory environment. Such contamination complicates the quantitative measurements of the initial amounts of phenolic compounds in foods, especially when they are present at low nanogram levels.

Methods

The method is an UPLC-ESI-MS/MS based multi-analyte method that applies multiple reaction monitoring (MRM) and isotope labeled internal standards for quantitative determination of phenolic compounds.

Results

This poster focuses on analysis of a selected number of most relevant phenolic compounds that can be found in foods, as well as approaches to minimize and control the external contamination from the laboratory environment during sample handling and analysis. The method is not validated yet and the findings from pre-validation measurements of phenols in some food samples, after preventive measures to reduce the external contamination sources in our laboratory were implemented, are presented. The levels of the found compounds were clearly higher than in the reagent blank.

Conclusions

According to the present experience, the analysis of phenolic EDCs at low concentration levels in mixed and complex food matrices using regular analytical methods and laboratory equipment seems to be practically unfeasible until all the environmental sources of these chemicals are identified and maximally reduced to minimize their recurrence and impact on the analysis in order to ensure the validity of the data.

Novel Aspect

The method used to analyze phenolic compounds in food is fast and straightforward. It allows analyzing different groups of phenolic compounds in the same analysis. The method also comprises important measures that have to be considered to minimize external contamination during the sample pre-treatment and thereby ensure contamination free extraction and reliability of the measurements.

References

1. Basile T, Petrella A, Petrella M, et al. Review of Endocrine-Disrupting-Compound Removal Technologies in Water and Wastewater Treatment Plants: An EU Perspective. *Industrial and Engineering Chemistry Research*, 50, 8389-8401, 2011.
2. Calafat A, Wong L-Y, Ye X, et al. Concentrations of the Sunscreen Agent Benzophenone-3 in Residents of the United States: National Health and Nutrition Examination Survey 2003–2004. *Environmental Health*

Perspectives, 116, 893-897, 2008.

3. Adamusova H, Bosakova Z, Coufal P, et al. Analysis of estrogens and estrogen mimics in edible matrices – A review. *Journal of Separation Science*, 37, 885-905, 2014.

WP-239 / A STUDY ON MIGRATION OF HEAVY METALS FROM CERAMICS AND GLASSWARE USING INDUCTIVELY COUPLED PLASMA MASS SPECTROMETRY(ICP-MS)

Kyung Chul Cho - Y. E. Cho - Yong Chjun Park - Dong me Choi

Gyeongin Regional Food and Drug Administration, Center for Food and Drug Analysis, Incheon, Republic of Korea

Keywords: ICP-MS, heavy metal, glassware, ceramics

Introduction

When ceramics and glassware come in contact with food, a migration of heavy metals such as lead and cadmium may occur [1]. The presence of heavy metals in foodstuff represents a crucial problem around the world. The risk associated with the exposure to heavy metals in food products had brought widespread concern in human health. This study investigated the migration level of heavy metals (lead, cadmium, arsenic and etc.) from food contact articles into a food simulant.

Methods

All the food contact articles (80 ceramics and 80 glassware) were purchased in domestic markets around South Korea. The articles were eluted with a food simulant (4% acetic acid) for 24 hr at 25°C in laboratory oven [2,3]. Heavy metals including lead and cadmium that migrated from articles into a food stimulant were analyzed by Inductively Coupled Plasma Mass Spectrometry (ICP-MS).

Results

The migration levels of heavy metals in all the samples were compared with the migration limits of Ministry of Food and Drug Safety (MFDS) in South Korea.

Conclusions

Analytical methods were successfully validated by evaluating linearity, LOD, LOQ, precision (intra/inter-day) and recovery.

Novel Aspect

The result of this study can be used as valuable data for the safety control of the food contact materials.

References

1. Gonzalez-Soto, E., Gonzalez-Rodriguez, V., Lopez-Suarez, C., Castro-Romero, J. M., Perez-Iglesias, J., Fernandez-Solis, J. M., *Bulletin of Environmental Contamination and Toxicology*, 65, 598–603 (2000).
2. Food Contact Materials - Regulation (EC) 1935/2004
3. C. Simoneau, B. Raffael, S. Garbin, E. Hoekstra, A. Mieth, J. A Lopes, V. Reina, JRC Science for policy report (2016).

WPS-S10 - Atomic mass spectrometry: metallomics and isotope ratio mass spectrometry

Chairs: Ryszard Lobinski, Federica Camin

WP-240 / ISOTOPIC NICHE SPACE OF SONGBIRDS MIGRATING THROUGH THE SOUTHERN ALPS

Alessandro Franzoi⁽¹⁾ - *Federica Camin*⁽²⁾ - *Paolo Pedrini*⁽¹⁾ - *Luana Bontempo*⁽²⁾

MUSE, Sezione Vertebrati, Trento (1) - Department of Food Quality and Nutrition, Research and Innovation Centre, Fondazione Edmund Mach, Via E. Mach 1, 38010 San Michele all'Adige, Italy (2)

Keywords: stable isotope ratios, Passerines, Alps, isotopic niche

Introduction

Stable isotopes are a practical tool used to define animal movements and to delineate trophic structures in animal communities [1, 2]. A crucial aspect of spatial ecology is the definition of the degree of isotopic variability that can be measured locally in animal tissues, in order to carry out the correct probabilistic assignments to a geographic location, calibrating the isotopic discrimination between stable isotopes in trophic sources and animal tissues [3].

Methods

Stable isotope ratios of bio-elements (²H, ¹⁸O, ¹³C, ¹⁵N, ³⁴S) were determined in feathers of Passerines from the Italian Alps. The feathers were sampled on an altitudinal gradient in a restricted region from juvenile and post-breeding generations. Stable isotope ratios were determined by Isotope Ratio Mass Spectrometry after preparation of the feathers according to Bontempo et al. (2014) [4].

Results

Isotopic variations can be mainly related to species feeding habits [5], seasonal effect of regional climatic conditions [6], geographic and topographic patterns [7, 8], natural or anthropogenic inputs which perturbate the local isoscapes [9, 10]. Through this study we assessed the probable variability sources, taking into account the distinct stable isotope ratios separately. In particular, the local variability in the sampling sites within a specific area of Italian Alps was measured, comparing the isotopic composition in feathers of several common Passerine species sampled from different moult generations in moulting and fattening sites within or close to the breeding locations. Using the derived information it was possible to determine the specific regional multi-isotopic variability, as well as to define seasonal isotopic variability related to trophic and habitat niches.

Conclusions

The local variability of the isotopic composition within the sampling sites, comparing different species and different feathers generations was measured. In particular, the regional multi-isotopic variability for Passerines according to trophic sources and animal tissues was defined. This information is crucial in depicting the isotopic niche in local birds of Central Alps, deepening the knowledge about isotopic variability related to topographic, environmental and seasonal perturbations.

Novel Aspect

Isotopic variability in avian keratinous tissues was estimated with a multi-isotopic approach in a narrow Alpine region. The concept of isotopic niche was deepened.

References

1. Hobson K. A. *Divers Distrib*, 11,157–164 (2005).
2. Layman C. A., Arrington D. A., Montana C. G., Post D. M., *Ecology*, 88,42–48 (2007).
3. Wunder M. B., Norris D. R., In: Hobson K. A., Wassenaar L. I. (eds) *Tracking animal migration with stable isotopes*. Elsevier, 107–128 (2008).
4. Bontempo L., Ceppa A., Ziller L., Pedrini P., Hobson K. A., Wassenaar L. I., Camin F., *methods Ecol Evol*, 5, 363-371(2014).
5. Hobson K. A., Van Wilgenburg S. L., Faaborg J., et al, *J F Ornithol*,85,237–257 (2014).

6. Bortolotti L. R., Clark R. G., Wassenaar L. I., *Ecol Appl*,23,110–121 (2013).
7. Bowen G. J., Wassenaar L. I., Hobson K. A., *Oecologia*, 143,337–348 (2005).
8. Newsome S. D., Sabat P., Wolf N., et al. *Ecosphere*, 6,1-18 (2015)
9. Hebert C. E., Wassenaar L. I., *Wildl Soc Bull*, 33,92–102 (2005).
10. Yi-Balan S. A., Amundson R., Buss H. L, *Geochim Cosmochim Acta*, 142,411–428 (2014).

WP-241 / EXPERIMENTAL EVALUATION OF TIMS-BASED DETECTION METHODS FOR ISOTOPIC ANALYSIS OF URANIUM AT ULTRA-TRACE LEVEL

Jong-Ho Park⁽¹⁾ - ***Kahee Jeong***⁽²⁾

Korea Atomic Energy Research Institute, Nuclear Chemistry Research Division, Daejeon (1) - Central Research Institute of Korea Hydro & Nuclear Power, Chemistry & Environmental Group, Daejeon (2)

Keywords: TIMS, Isotopic Analysis, Uranium, Detection Method, Nuclear Safeguards

Introduction

Thermal ionization mass spectrometry (TIMS) is one of the most accurate and precise techniques for isotopic analysis, and suitable for monitoring undeclared nuclear activities[1-3]. The analytical performance of the three TIMS-based detection methods was experimentally evaluated for isotopic analysis of uranium at ultra-trace levels, in terms of analytical accuracy, precision, and measurement uncertainty.

Methods

Isotopic measurement of uranium at 1 ng, 100 pg, 30 pg, 5 pg, and 1 pg levels were performed by TIMS (TRITON Plus, Thermo) using three detection methods; Dynamic, multi-dynamic, and static methods. Analytical performance, such as accuracy, precision, and measurement uncertainty, was experimentally evaluated depending on the detection methods.

Results

The TIMS analysis for the ten replicated U samples with amounts of 1 ng and 100 pg showed no significant improvement in analytical performance irrespective of the detection method adopted, whereas slight improvement was observed for the analysis of 30 pg uranium using the multi-dynamic and static detection methods. Considering the analytical performance and the easiness of detector calibration, the preferred detection methods for 1 ng and 100 pg U is the multi-dynamic method, while that for 30 pg U is the dynamic method. The static detection method offers the greatest accuracy and precision, and the smallest uncertainty for TIMS measurements of 5 pg and 1 pg of U due to the greater detection sensitivity of ion counters than faraday cups, the elimination of ion signal drift, and the large number of valid data sets in a measurement.

Conclusions

This study deals with the evaluation of three detection methods (multi-dynamic, dynamic, and static) in the isotopic analysis of ultra-trace amounts of uranium using thermal ionization mass spectrometry. For samples containing 5 pg and 1 pg uranium, the static detection method provided the best results in terms of accuracy, precision, and measurement uncertainty.

Novel Aspect

This study shows that the static detection method of TIMS can be applied to particle analysis of environmental samples for nuclear safeguards.

References

1. Donohue D., *Anal. Chem.*, 74, 28A (2002)

2. Heumann K., Eisenhut S., Gallus S., Hebeda E., Nusko R., Vengosh A., Walczyk T., Analyst, 120, 1291 (1995)
3. Stetzer O., Betti M., Geel J., Erdmann N., Kratz J., Schenkel R., Trautmann N., Nucl. Inst. Meth. Phys. Res. A, 525, 582 (2004)

WP-242 / C AND H STABLE ISOTOPE RATIO ANALYSIS USING GC-IRMS FOR VANILLIN AUTHENTICATION

Matteo Perini⁽¹⁾ - **Silvia Pianezze**⁽¹⁾ - **Federica Camin**⁽²⁾

Technology Transfer Centre, Fondazione Edmund Mach, via E. Mach 1, 38010 San Michele all'Adige (TN), Italy (1) - Fondazione Mach, Department of Food Quality and Nutrition, Research and Innovation Centre, S. Michele all'Adige (TN), Italy (2)

Keywords: stable isotope ratios, vanillin, authentication, gas chromatography

Introduction

Vanilla extracts are widely used as flavouring ingredients in foods and beverages and as aromatic compounds in perfumes and pharmaceuticals. Due to the high production cost of high-quality natural extracts from *Vanilla planifolia*, synthetic or natural identical biosynthetic vanillin (which derives from natural precursors such as guaiacol, ferulic acid, eugenol and lignin) are often used as a substitute for authentic natural vanillin [1].

Methods

We combined analysis of $^{13}\text{C}/^{12}\text{C}$ (expressed as $\delta^{13}\text{C}$) values and $^2\text{H}/^1\text{H}$ (expressed as $\delta^2\text{H}$) ones using gas chromatography coupled to isotope ratio mass spectrometry (GC-IRMS) [2]. 16 authentic samples of *Vanilla planifolia*, 16 natural identical, 5 synthetic vanillin and 20 commercial extracts were considered.

Results

The carbon isotope ratio values allow discriminating between natural and natural identical or synthetic vanillin (averagely -19‰ vs -30‰). Authentic natural vanillin from *Vanilla planifolia* and natural identical vanillin are characterised by $\delta^2\text{H}$ values much lower than synthetic vanillinones (averagely -50‰ vs 73‰). The isotopic values of all the commercial extracts declared to be from *Vanilla planifolia* (N=20), had $\delta^{13}\text{C}$ values within the typical range of natural vanillin, but $\delta^2\text{H}$ values outside the range and more similar to that of synthetic vanillin. Our results show how the stable isotope ratio analysis (SIRA) of $^{13}\text{C}/^{12}\text{C}$ is no longer sufficient to discover vanillin adulteration, due to the practice of adding ^{13}C to the methyl site of synthetic vanillin [3], and the new analysis of $\delta^2\text{H}$ is mandatory.

Conclusions

The combination of $\delta^{13}\text{C}$ with $\delta^2\text{H}$ GC-IRMS analysis of vanillin can therefore be proposed as a suitable tool to improve the detection of vanilla extract adulteration.

Novel Aspect

Stable isotope analysis of $^2\text{H}/^1\text{H}$ can be a rapid and effective method to identify the natural or synthetic origin of vanillin.

References

1. Bensaid FF, Wietzerbin K, Martin GJ. Authentication of natural vanilla flavorings: isotopic characterization using degradation of vanillin into guaiacol. J Agric Food Chem. 2002;50: 6271–6275.
2. Greule M, Tumino LD, Kronewald T, Hener U, Schleucher J, Mosandl A, et al. Improved rapid authentication of vanillin using $\delta^{13}\text{C}$ and $\delta^2\text{H}$ values. Eur Food Res Technol. 2010;231: 933–941.
3. Krueger DA, Krueger HW. Detection of fraudulent vanillin labeled with carbon-13 in the carbonyl carbon. J Agric Food Chem. 1985;33: 323–325.

WP-243 / STABLE ISOTOPE RATIO ANALYSIS FOR AUTHENTICATION OF CHITOSAN

Matteo Perini - Tiziana Nardin - Matteo Venturelli - Silvia Pianezze - Roberto Larcher

Technology Transfer Centre, Fondazione Edmund Mach, via E. Mach 1, 38010 San Michele all'Adige (TN), Italy

Keywords: stable isotope ratios, chitosan, authentication, animal or fungi origin

Introduction

Chitosan is a linear polysaccharide with a number of possible uses (e.g. in medicine). It is produced by deacetylation of chitin, which is the structural element of exoskeleton of crustaceans. Due to the presence of allergens in seafood, its production from fungi has gained increasing attention. To identify its origin, residual glucans, viscosity and settled density are used, but often not effective. This calls to the development of new methods, such as the stable isotope ratio analysis.

Methods

Following the prescription of OIV International Oenological Codex [1], residual glucans, viscosity and settled density were measured to identify the crustacean or fungi origin of 20 chitosan samples. Moreover, the stable isotope ratios (SIRs) of four bioelements - hydrogen, carbon, nitrogen, oxygen - were analysed in the raw product using Isotope Ratio Mass Spectrometry.

Results

The biosynthetic pathway of the chitosan's precursor chitin affects its isotopic composition. The samples identified by the OIV official methods as from crustaceans showed lower $\delta^{13}\text{C}$ (averagely -20.8‰) and higher $\delta^{15}\text{N}$ values (-2‰). This fits with the marine origin of the font (the exoskeleton of crustaceans) characterized by low $\delta^{13}\text{C}$ value (around -20‰) and high positive $\delta^{15}\text{N}$ value (around +7‰)[2]. The low $\delta^{15}\text{N}$ values here found are probably due to the addition to the yeast broth of nitrogen synthetic sources such as urea with low $\delta^{15}\text{N}$ values [3]. The samples identified of fungi origin based on the OIV prescriptions, are characterized by higher $\delta^{13}\text{C}$ (-13.8‰) and lower $\delta^{15}\text{N}$ values (-4.1‰). The higher $\delta^{13}\text{C}$ values of these chitosan samples are explained on the basis of the fact that sugar cane juice and molasses with typical $\delta^{13}\text{C}$ around -11‰ [4] are used as economic substrates in fungi chitosan production [3]. D/H and $\delta^{18}\text{O}$ are lower in chitosan from fungi probably due to the lower value of the fermentation water.

Conclusions

This work considers a robust and effective method, based on the determination of stable isotope ratios of bioelement, that can be used to identify the origin of chitosan.

Novel Aspect

Stable isotope analysis could be a rapid method to identify the crustaceans or fungi origin of chitosan.

References

1. Resolution OENO 368/2009, Chitosan, OIV International Oenological Codex
2. Boutton T. W., Carbon Isotope Techniques. Academic Press, 173 (1991).
3. Kurtzman C. P., Fell J. W., Boekhout T., Robert V. methods for isolation, phenotypic characterization and maintenance of yeasts. Elsevier, 1, 87 (2011)
4. Bauer-Christoph C., Wachter H., Christoph N., Rossmann A., Adam L. Z Lebensm Unters Forsch A, 204, 445 (1997).

WP-244 / DEVELOPMENT OF NOVEL COMPACT CARBON DATING INSTRUMENT BASED ON COLLINEAR RESONANCE IONISATION SPECTROSCOPY (CRIS)

Giles Edwards - Kieran Flanagan - Ben Cooper - Sultan Alsufyani - Christopher Ricketts - Cory Binnersley

The University of Manchester, School of Physics and Astronomy, Manchester, United Kingdom

Traditional carbon dating techniques based on isotope ratio accelerator mass spectrometry are expensive, experimentally complex and access to instrumentation is often associated with excessive lead times. The use of CRIS as an alternative method affords unparalleled interference suppression when compared to RIMS for example. This is achieved by effectively decoupling the inlet electron cyclotron resonance ion source from the resonance ionisation spectroscopy by means of a charge exchange cell. The instrument employs a variety of ion optics including initial multi-collector magnetic sector mass selection and a linear ion trap to increase duty cycle; the sample inlet is a flash GC nitrogen/carbon elemental analyser. The CRIS technique was initially developed at CERN in the ISOLDE facility for the analysis of exotic nuclei. The instrument at the University of Manchester is configured for the analysis of carbon isotopes but can easily be tuned for other nuclei and offers a state-of-the-art alternative technique for various forensic, dating and food security applications.

THURSDAY, August 30th, 2018

ThPS-S01 - Glycomics, saccharides and glycoconjugates

Chairs: Catherine Costello, Joseph Banoub

ThP-1 / REDUCTIVE ALKALINE RELEASE OF N-GLYCANS GENERATES A VARIETY OF UNEXPECTED, USEFUL PRODUCTS

Rudolf Figl - Friedrich Altmann

University of Natural Resources and Life Sciences, Chemistry / Glycobiology, Vienna

Keywords: N-glycans, alkaline release, PGC-LC-ESI-MS

Introduction

Release of O-glycans by reductive β -elimination has become routine in many glyco-analytical laboratories and concomitant release of N-glycans has repeatedly been observed. Revisiting this somewhat forgotten mode of N-glycan release revealed that all kinds of N-glycans including oligomannosidic and complex-type N-glycans from plants with 3-linked fucose and from mammals with or without 6-linked fucose and with sialic acid could be recovered.

Methods

Typical conditions for the alkaline borohydride treatment of glycopeptides were varied in case of temperature and concentration of chemicals. Following PGC SPE cleanup, yield and chemical nature of the products were investigated using MALDI-TOF MS and PGC-LC-ESI-MS. Reducing glycans were labelled with 2-aminobenzamide, reducing glycosylamine and 1-amino-1-deoxy-glycitol with InstantPC™ dye. The fluorescent products were separated by HILIC-FLD-HPLC.

Results

Even after 16 h incubation in 1 M sodium borohydride, a large part of the glycans occurred in reducing form. Moreover, about one third emerged in the form of the stable amino-functionalized 1-amino-1-deoxy-glycitol. Relative yields of the products reducing glycosylamine, reducing N-glycan, 1-amino-1-deoxy-glycitol or glycitol could be controlled by the release conditions.

Conclusions

Chemical release of N-glycans constitutes a cost-saving alternative to enzymatic hydrolysis for the preparation of precursors for the production of reference compounds for various formats of N-glycan analysis. Moreover, it allows to obtain a stable amino-functionalized glycan derivative, which can be

employed to construct glycan arrays or affinity matrices.

Novel Aspect

While the principle observation of the release of *N*-glycans by alkali under reducing conditions was made decades ago, the very nature of the generated products has not been recognized so far.

References

1. R. Figl, F. Altmann, *Proteomics* 2018, 0, 1700330.

ThP-2 / ANALYSIS OF MOUSE BRAIN *N*-GLYCANS BY (PGC-)LC-ESI-MS USING ISOTOPICALLY LABELED STANDARDS

Clemens Grünwald-Gruber - Andreas Thader - Jonathan Urteil - Daniel Maresch - Friedrich Altmann

University of Natural Resources and Life Sciences, Department of Chemistry, Vienna

Introduction

N-glycosylation is the most complex and heterogeneous post translational modification. It is important for protein folding, stability and biological interactions of many proteins. The large number of possible isomers yielding from one glycan composition necessitates powerful separation techniques (porous graphitic carbon LC) and the use of defined standards combined with MS/MS.

Methods

Proteins of mouse brain tissue was purified, treated with PNGase F to release the *N*-linked glycans and reduced with sodium borohydride. The glycans were purified and directly subjected to a PGC-LC-ESI-MS system (QTOF instrument). Prior to the measurement defined isotopically labeled *N*-glycan standards were added to identify and quantify the isomers occurring. The incorporation of different numbers of heavy isotopes at varying positions ($^{13}\text{C}_6$ - and $^{13}\text{C}_1$ galactose, reduction with sodium borodeuteride or deuterio acetyl groups of GlcNAc) allowed the assembly of internal multiplex standards mixtures.

Results

The major glycan was found to be a bisecting, fucosylated GnGn (GnGnF (bi) -proglycan nomenclature, www.proglycan.com). All main isoforms could be identified and relatively quantified. The ability of PGC to discriminate about 30 different isobaric mono-fucosylated *N*-glycans was shown (H5N4F1). Absolute quantification could be furthermore achieved by the quantitation of one ^{13}C labelled standard by amino sugar analysis. This one internal standard allowed for the calculation of the other glycoforms by the use of instrument dependent correction factors (calculated ratios known from an equimolar standard mix [1])

Conclusions

A powerful method for the analysis of *N*-glycans in complex tissue samples is presented. The use of different labeled standards allowed identification and quantification of glycan isoforms in a single MS-run.

References

1. Grünwald-Gruber C, Thader A, Maresch D, Dalik T, Altmann F. Determination of true ratios of different *N*-glycan structures in electrospray ionization mass spectrometry. *Anal Bioanal Chem.* 2017 Apr 409, 2519-2530

ThP-3 / ASSIGNMENT OF LINKAGE, ANOMERIC CONFIGURATION, AND BRANCH LOCATION OF OLIGOSACCHARIDES VIA A LOGICALLY DERIVED MS^n APPROACH

***Hsu Chen Hsu*⁽¹⁾ - *Chia Yen Liew*⁽¹⁾ - *Shih-Pei Huang*^(1,2) - *Shang-Ting Tsai*⁽¹⁾ - *Chi-Kung Ni*^(1,3)**

Institute of Atomic and Molecular Sciences, Academia Sinica, Taipei (1) - Department of Chemistry, National Taiwan Normal University, Taipei (2) - Department of Chemistry, National Tsing Hua University, Taipei (3)

Keywords: carbohydrate, branched and linear oligosaccharides, structural determination

Introduction

Carbohydrates (or glycan) play influential roles in cellular recognition processes and many other biological processes. However, the structural determination of oligosaccharides remains an obstacle and constitutes a bottleneck in carbohydrate research. An approach for de novo structural determination of linear and branched oligosaccharides based on logically derived procedure of tandem mass spectrometry was demonstrate.

Methods

A commercial linear ion trap mass spectrometry (LTQ XL, Thermo Scientific) coupled with electrospray ionization was employed in this study. The selection of structural decisive fragment ions for the subsequent CID was guided using a procedure that we built from the understanding of the saccharide dissociation mechanism. The complete structural identification is obtained by combining the structural information from each step of CID.

Results

A specially prepared disaccharide CID spectrum database consisted of 16 mannobioses was build. Low-energy collision-induced dissociation (CID) of sodium ion adducts was used to facilitate the cleavage of desired chemical bonds during dissociation. The linkages, anomeric configurations, and branch locations of oligosaccharides are determined by comparing the CID spectra of oligosaccharide with the fragmentation patterns based on dissociation mechanism and our specially prepared disaccharide CID spectrum database. The applications of this method to structural determination of mannose-oligosaccharides were demonstrated. The structure of a given trisaccharides from 88 isomers and tetrasaccharide from 928 isomers can be determined by only 3-5 properly selected CID mass spectra according to the logical procedure we derived. The results show that this method provides a rapid structural identification procedure for mannose oligosaccharides and can be extended to glucose and galactose oligosaccharides as well.

Conclusions

We demonstrated a simple method that an oligosaccharde can be identified from hundreds isomers by only using 3-5 properly selected CID spectra. The required CID spectra for structural identification can be determined from the schemes we derived according to the dissociation mechanisms. The entire procedure for these four CID spectrum measurements only take few minutes, demonstrating this approach is a simple and rapid method for the structural determination of oligosaccharides.

Novel Aspect

The approach can be built in a computer controlled mass spectrometer and ultimately developed into a automatic procedure for structural identification of oligosaccharides.

References

1. Hsu H. C., Liew C. Y., Huang S. P., Tsai S. T., Ni C. K., J. Am. Soc. Mass Spectrom. 29, 470, (2018).
2. Tsai S. T., Chen J. L., Ni C. K., Rapid Commun. Mass Spectrom. 31, 1835, (2017).
3. Hsu H. C., Liew C. Y., Huang S. P., Tsai S. T., Ni C. K., Sci. Rep., 8,5562, (2018).

ThP-4 / FRAGMENTATION OF OLIGOSACCHARIDES AND SUPPRESSION OF METAL-SALT INDUCED ADDUCTS USING ELECTROSPRAY IONIZATION

Volker Iwan - Jürgen Grotemeyer

Keywords: Mass spectrometry, FT-ICR, Oligosaccharides, CID, alkaline adduction

Introduction

It is known that alkaline adduct-ions influence CID and PD induced fragmentation mechanisms of oligosaccharides [1,2]. For proteins the use of ammonium acetate to counteract cation induced signal suppression in electrospray ionization is published [3]. To increase the intensity of the non-adducted molecular ion the application of ammonium acetate to oligosaccharide containing samples is of interest.

Methods

The experiments were performed on an APEX Qe FT-ICR mass spectrometer, with a 9.4 T superconducting magnet (Bruker Daltonik, Bremen, Germany). Maltopentaose (MP) was dissolved in MeOH/H₂O and ionized using an Apollo III ESI-source. To analyse the suppression of alkali-adducts NH₄Ac was added to the samples. Fragmentation of the molecular cations were achieved by SORI-CID with argon as collision gas.

Results

Several spectra of maltopentaose with different concentrations of NH₄Ac could be recorded. The S/N-ratio of the non-adducted molecular ion could be increased by nearly sixty times compared to the measurement without the use of ammonium acetate. However, no complete suppression of the adducted molecules could be achieved. The use of NH₄Ac leads to an increased fragmentation - similar to CID experiments. Nevertheless, it allowed the comparison of SORI-CID indicated fragmentation behaviour of the different species. The spectra of the protonated species show no cross-ring cleavage, but loss of H₂O and glucose. Compared to that the CID spectra of the Na-adducts show cross-ring cleavages, also in combination with loss of glucose. In this regard the study of the effect on chromophore labelled maltopentaose and subsequent CID/PD fragmentation would be part of interest and is indeed part of ongoing research.

Conclusions

Using ammonium acetate to reduce alkaline adduct ions of oligosaccharides was successfully applied to the example of MP, whereby MS² experiments of the [M+H]⁺ species could be performed. Doing this, no cross-ring cleavages are observable. Due to the symmetry of MP, an assignment of the fragments to the reducing or non-reducing end is not possible. Because of the precipitation of acetates and ammonium salts, this method requires a higher cleaning effort.

Novel Aspect

The use of ammonium acetate for the suppression of alkaline adduct-ions to improve the S/N-ratio of the [M+H]⁺ species measuring oligosaccharides.

References

1. Hahn A., Grotemeyer J., Anal. Methods, 5, 503-507 (2013).
2. Cancilla M., Wong A., Voss L., Lebrilla C., Analytical Chemistry, 71, 3206 (1999).
3. Iavarone A., Williams E., Analytical Chemistry, 76, 3944 (2004).

ThP-5 / MALDI-MS ANALYSIS OF SUCROSE USING A CHARCOAL MATRIX

Jeongkwon Kim

Chungnam National University, 99 Daehak-Ro, Yuseong-Gu, Daejeon

Keywords: sucrose, charcoal, MALDI, glycosidic bond, cationization agent

Introduction

Matrix-assisted laser desorption/ionization mass spectrometry (MALDI-MS) analysis of low molar mass

compounds using a conventional matrix, such as 2,5-dihydroxybenzoic acid (DHB) or α -cyano-4-hydroxycinnamic acid (CHCA) is not commonly performed due to matrix-related peaks [1, 2]. For the effective MALDI-MS analysis of sucrose whose average molar mass is only 342.3 Da, a charcoal matrix with different cationization agents was utilized.

Methods

Different cationization agents such as Li^+ , Na^+ , K^+ , Rb^+ , Ag^+ , and Cs^+ were used for the MALDI-MS analysis of sucrose with a charcoal matrix. The matrix (1 μL) was placed on a MALDI plate and dried at room temperature. Then, 1 μL of 100 mM sucrose stock solution was deposited on top of the matrix. To compare the effectiveness of cationization agents, each additive cation (1 μL) was loaded on top of the dried spots. Mass spectra were obtained in the positive ion reflection mode.

Results

A higher cation-adducted sucrose peak with significantly reduced interference was observed with a charcoal matrix rather than DHB or CHCA. However, the charcoal matrix caused glycosidic bond cleavage, resulting in sucrose fragment peaks. These could not be removed, even when we reduced the laser intensity. The degree of sucrose fragmentation was inversely related to the size of the cation additive. More sucrose fragmentation occurred with small cationization agents (Li^+ , Na^+ , or Ag^+), while little fragmentation occurred with relatively large cationization agents (K^+ , Rb^+ , or Cs^+).

Conclusions

Charcoal has a higher energy transfer efficiency than DHB or CHCA. This may explain the increase in sucrose peaks and fragmented peaks observed when sucrose was analyzed with a charcoal matrix. The degree of sucrose fragmentation was inversely related to the size of the cation additive. The highest degree of sucrose fragmentation was observed when the cation additive was Li^+ , followed by Ag^+ , Na^+ , and K^+ . No fragmentation was observed with additives Rb^+ and Cs^+ .

Novel Aspect

Charcoal matrix was utilized for the analysis of sucrose in MALDI-MS with different cationization agents such as Li^+ , Na^+ , K^+ , Rb^+ , Ag^+ , and Cs^+ .

References

1. Kim J., Mass Spectrom. Lett. 6, 27 (2015).
2. Kim J. -I., Park J. -M., Noh J. -Y., Hwang S. -J., Kang M. -J., Pyun J. -C., Chemosphere 143, 64 (2016)

ThP-6 / MULTIPLEXED ISOTOPE LABELED STANDARDS FOR ISOMERIC ANALYSIS OF N-GLYCANS

Daniel Maresch - Clemens Grünwald-Gruber - Andreas Thader - Friedrich Altmann

University of Natural Resources and Life Sciences, Department of Chemistry, Vienna

Introduction

Porous graphitic carbon (PGC) columns have proven extraordinarily useful for the LC-ESI-MS analysis of *N*- and *O*-glycans. Their shape selectivity has been demonstrated for different sialic acid and galactose linkages. The ability of PGC to discriminate about 35 different isobaric mono-fucosylated *N*-glycans will be shown. Application to mouse brain *N*-glycans revealed them to contain hybrid type glycans with fucose at the core and Le x (but not Le a or blood group H) determinants.

Methods

The notoriously poor reproducibility of retention times on PGC prompted us to generate internal standards with different mass increments due to thoughtful incorporation of one or two molecules of either $^{13}\text{C}_6$ - or $^{13}\text{C}_1$ -galactose, of BH_4 or BD_4 reduction and finally de-*N*-acetylation and re-*N*-acetylation using stable-isotope labeled acetic anhydride. The obtained multiplexed reference glycans can easily be told apart by

their specific mass increment.

Results

Since expression of retention times as glucose units is not useful for LC on PGC, we propose to scale any chromatogram to a reference run and express the retention in the form of virtual minutes (vimin). Absolute quantification was furthermore achieved by the quantitation of one $^{13}\text{C}_6$ labelled standard by amino sugar analysis. This one internal standard allowed for the calculation of the other glycoforms by the use of instrument dependent correction factors (calculated ratios known from an equimolar standard mix [1]).

Conclusions

A strategy is proposed by which a set of multiplexed standards can serve to assign the structure of any *N*-glycan structure already by retention time, thus providing a new level of confidence in the analysis of structural isomers that supports and expands the conclusions drawn from MS/MS spectra.

References

1. Grünwald-Gruber C, Thader A, Maresch D, Dalik T, Altmann F., Anal Bioanal Chem. 409, 2519-2530 (2017)

ThP-7 / *N*-GLYCAN ANALYSIS OF VARIOUS CHLORELLA STRAINS BY MALDI-TOF AND LC-ESI-MS

Reka Tunde Mocsai - Rudolf Figl - Elisabeth Svehla - Friedrich Altmann

University of Natural Resources and Life Sciences, Department of Chemistry, Vienna

Keywords: Chlorella, MALDI-TOF, LC-ESI-MS, *N*-glycosylation

Introduction

Microalgae of the genus *Chlorella* are candidates for the production of polyunsaturated fatty acids, various dietary supplements and they are also debated as novel expression systems for recombinant proteins. Our work shows that different *Chlorella* species display widely varying *N*-glycan profiles with unique masses not known from embryophytes (land plants).

Materials and Methods

Total *N*-glycan spectra of a panel of *Chlorella vulgaris* and *Chlorella sorokiniana* strains as well as of a broad selection of commercial *Chlorella* tablets - many of them labeled with the outdated taxonomical name *Chlorella pyrenoidosa* - were investigated by MALDI-TOF-MS and LC-ESI-MS.

Results

Most samples typed as *C. vulgaris* contained oligomannosidic glycans with several *O*-methyl groups, while the *C. sorokiniana* culture strains exhibited two differing profiles with pentose containing glycans. The commercial samples -almost irrespective of their label- yet again occurred as four distinct structural groups. One resembled *C. vulgaris*, two showed rather complex patterns of methylated glycans, and one contained just one very dominant structure that turned out to contain three arabinoses, one of them methylated and a galactofuranose linked to the second *N*-acetylglucosamine residue.

Conclusions

Chlorella N-glycans differ greatly from plant and vertebrate *N*-glycans by the lack of fucosylation, the presence of arabinose and galactose residues and the prominent *O*-methylation of sugars. These properties foresee a large variety of unknown glycosyltransferases and potentially immunogenic glycan structures. The occurrence of groups with clearly distinct *N*-glycan patterns suggests that *N*-glycan profiles could constitute a valuable tool for taxonomic classification as well as for determination of the origin of *Chlorella* food supplements.

ThP-8 / SEMI-TARGETED ANALYSIS OF METABOLITES IN SAPONARIA OFFICINALIS BY LC-MS

Katarzyna Pawlak⁽¹⁾ - ***Katarzyna Lech***⁽¹⁾ - ***Kamil Wojciechowski***^(1,2) - ***Ilona Goral***⁽³⁾ - ***Piotr Bartosz***⁽¹⁾ - ***Piotr Kosior***⁽¹⁾ - ***Tomotheus Kiehl***⁽¹⁾ - ***Anna Ornoch***⁽¹⁾ - ***Kinga Kaczmarczyk***⁽¹⁾

Warsaw University of Technology, Faculty of Chemistry, Chair of Analytical Chemistry, Warsaw (1) - SapoLabs Ltd, Chair of Medical Biotechnology, Warsaw (2) - Warsaw University of Technology, Faculty of Chemistry, Chair of Medical Biotechnology, Warsaw (3)

Keywords: semi-targeted, metabolomics, saponins, plants, *S. officinalis*

Introduction

Saponaria officinalis is a common perennial plant from the Caryophyllaceae family. It is well-known for its detergent property and has been used as a soap in ancient times, later also as an expectorant in bronchitis and in rheumatic disorders. Decocts and infusions of this plant can be used as an emulsifier for food and cosmetic products (masks, toothpastes, shampoos). However, the fact that its properties depend on the organ and time of harvest is often overlooked.

Methods

Saponins were tracked by LC-MS system using semi-targeted metabolomic approach. Changes in the metabolome compositions were registered by following 400 transitions in negative ion mode by 6460 ESI MS/MS system coupled to RPLC. For each parent ion (selected from the plant's analysis in the scanning mode, range 350 – 1700 *m/z*) at least two product ions were selected to track changes of metabolomics profiles for different organs and harvest time.

Results

Different triterpenoid saponins, saponariosides A-M named and classified by Kazuo Koike and Zhonghua Jia [1], were detected in extracts from the whole plants of *Saponaria officinalis*. Throughout the ESI-MS profiling of RPLC eluate, it was recognized that lower molecular weight saponins of *S. officinalis* were presenting the very polar fractions eluted faster. Very high diversity of the structures was observed as different sugars were found to be attached to aglycone, such as pentose, hexose, and deoxy-sugars. This was indicated as the difference between *m/z* values for signals such as 132, 162 and 146 for deoxyhexosyl moieties. Cluster analysis of the collected peak areas for individual saponins were visualized as a dendrogram and heat map in respect to time of harvest. We found that flowering and formation of seeds influences significantly the relative amounts of different saponariosides. The observed changes in composition show good correlation with parameters illustrating the ability of extracts to form stable and dense foam.

Conclusions

The HPLC-MS/MS analysis allowed us to detect and identify saponins, and perform their semi-quantification. Developed LC-MS/MS platform allowed to track changes in saponin's composition for different organs obtained for *S. officinalis* harvested in the period from June to October.

Novel Aspect

For the first time, the stage of the development and organ-specificity of the crop were included in studies of *S. officinalis*.

References

1. Koike K., Jia Z., Nikaido T., J. Nat. Prod., 62, 1655 (1999)

This work was financially supported by NCBiR project POIR.01.02.00-00-0005/16

ThP-9 / SEMI-AUTOMATIC SITE SPECIFIC QUANTITATIVE ANALYSIS OF HIGH MANNOSE AND HYBRID TYPE GLYCOSYLATION IN LIVER TISSUE ASSOCIATED WITH HCC AND CIRRHOSIS

Miloslav Sanda - Julius Benicky - Zuzana Brnakova - Lihua Zhang - Radoslav Goldman

Georgetown University, Department of Oncology, Washington, DC

Keywords: Glycoproteomics, Oligomannoseglycans, Hepatocellular carcinoma, Glycopeptide quantification, Liver disease

Introduction

Glycosylation changes of serum glycoproteins associated with liver disease progression were described previously. Most of the studies were profiling detached glycans from crude human sera or isolated proteins from patient sera. Recent studies described more detailed site specific glycosylation analysis of isolated glycoproteins or site specific analysis of particular glycan modifications such as core-fucosylation after glycosidase treatment. Our study is focused on the site specific analysis of incompletely processed glycan structures as oligomannose and hybrid glycans in liver tissue associated with liver disease progression (cirrhosis, HCC). Lectin affinity enrichment combined with HR mass spectrometry analysis allowed us to identify large set of high mannose, hybrid and paucimannose glycopeptides in complex liver tissue glycoproteome.

Methods

Human Liver tissue biopsies divided to tumor and adjacent portions were analyzed as whole tissue lysate and lectin enriched glycopeptides. Tryptic digestion including reduction and alkylation steps was followed by ConA lectin enrichment of oligomannose and hybrid glycopeptides. Whole liver proteome and enriched glycopeptides were analyzed by HR LC-MS/MS using Orbitrap Fusion Lumos mass spectrometer working in DDA mode coupled with 3000 RSLC nano-chromatography (Dionex). Peptides and glycopeptides were separated on the 50 cm C18 capillary column using 90 min acetonitrile gradient. Data were processed using Proteome discoverer (Thermo), Byonic and GPS software. All glycopeptide spectra were confirmed manually.

Results

Complex site specific glycosylation analysis as one of the most difficult analytical approach was never completely solved because of glycopeptide mass spectrometric fragmentation similarity. Using combination of selective structure enrichment, HR mass spectrometric analysis, chromatography reproducibility we were able to identify large set of high mannose and hybrid glycopeptides (2429) in human liver tissue. Byonic, GPS software and retention behavior prediction help us to identify additional high mannose and hybrid glycopeptides which were masked in crude liver tissue analysis by oligomannosylated peptides from abundant complex type glycopeptides. Large portion of identified ConA enriched fraction were trimmed paucimannose glycostructures (734) related to increased lysosomal activity. Some of the paucimannose glycostructures of lysosomal glyco-proteins for example prosaposin, TPP, Cathepsins or LAMP are elevated in liver HCC.

Conclusions

We have identified and quantified 2429 high mannose and hybrid glycopeptides in 6 human liver biopsies samples divided to tumor and adjacent tissue. We have detected elevation of proteins and glycoforms in tumor related to lysosomal activity.

Novel Aspect

High mannose and hybrid type glycoproteomics related to liver disease.

ThP-10 / INTEGRATING ION SPECTROSCOPY WITH MASS SPECTROMETRY FOR THE STRUCTURAL ANALYSIS OF CARBOHYDRATES

Baptiste Schindler - Gina Renois-Predelus - Isabelle Compagnon

Institut Lumière Matière, Université de Lyon, Lyon

Keywords: structural analysis, IRMPD spectroscopy, hyphenated MS-based approach, carbohydrate sequencing

Introduction

Sequencing techniques have been established for proteins and DNA and have revolutionized modern biology but routine techniques do not exist for glycans. In particular, top down analysis of glycans by mass spectrometry is often ambiguous due to the presence of various isomerisms. This lack of structural tools is identified as a crucial bottleneck, limiting the full development of glycosciences.

Methods

Our instrument combines Mass Spectrometry and Vibrational Spectroscopy (InfraRed Multiple Photon Dissociation) in the 3 μm spectral range. It offers the structural resolution of spectroscopic techniques for mass-selected ions and operates at typical MS conditions, *i.e.* it requires significantly less sample than traditional spectroscopy.

Results

First, we have shown that the IR signature in the 3 μm range obtained with our instrument is a powerful metric to resolve simultaneously all isomerisms: the nature of the monosaccharide content, the regiochemistry of the glycosidic bond, the anomer stereochemistry (α or β), the structure of the monomeric ring (pyranose or furanose), and the position of modifications such as sulfation or acetylation. Then the conservation of the molecular structure of a precursor ion within MS fragments has been revealed on disaccharides, opening the way to top-down, *de novo* MS analysis. [1] Following this demonstration, we have established a set of carbohydrate sequencing rules using a combination of Mass Spectrometry and IRMPD Spectroscopy. Finally, we have applied our approach to the determination of the sequence of various classes of saccharides. [2]

Conclusions

Without chemical derivatization of the glycan and with a minimal database of monosaccharide standards, the monosaccharide content and the stereochemistry of any oligosaccharide can be retrieved using our hybrid analytical approach based on MS-sequencing strategies and the structural resolution of carbohydrate isomerisms by IRMPD spectroscopy.

Novel Aspect

Our approach is to combine the best of two worlds: Mass Spectrometry and Spectroscopy, for the identification of glycan isomers.

References

1. Schindler B., Barnes L., Renois G., Gray C., Chambert S., Fort S., Flitsch S., Loison C., Allouche A. -R., Compagnon I., Nature Communications, 8, 973 (2017).
2. Schindler B., Renois-Predelus G., Bagdadi N., Melizi S., Barnes L., Chambert S., Allouche A. -R., Compagnon I., Glycoconj J, 34, 421-425 (2017).

ThP-11 / ISOMER-SPECIFIC MRM PROFILING OF SERUM N-GLYCOME FOR BEHCET DISEASE MONITORING

***Nari Seo*⁽¹⁾ - *Kyoung Heon Kim*⁽²⁾ - *Joong Kyong Ahn*⁽³⁾ - *Hoon Suk Cha*⁽³⁾ - *Jaehan Kim*⁽⁴⁾ - *Hyun Joo An*⁽¹⁾**

Graduate School of Analytical Science and Technology, Chungnam National University, Daejeon (1) - Department of Biotechnology, Korea University, Seoul (2) - School of Medicine, Sungkyunkwan University, Seoul (3) - Department of Food and Nutrition, Chungnam National University, Daejeon (4)

Keywords: Behcet disease, *N*-Glycan, Isomer, Quantitation, MRM-MS

Introduction

Glycosylation is widely recognized as the powerful resource of biomarkers for cancers and infectious diseases. Behcet disease (BD), a typical immune disease, is a vasculitis disorder characterized by oral sores and arthritis. Until now several diagnostic criteria of BD have been proposed however, there is no reliable diagnosis/treatment method. Here, we developed glycomic-based novel approach for BD monitoring using PGC-MRM MS.

Methods

Sera of Behcet patients (n=47) and health controls (n=47) were obtained from Samsung Medical Center in Korea. Each sample was enzymatically treated with PNGase F to release *N*-glycans. Native *N*-glycans were enriched by PGC-SPE prior to MS analysis. Targeted acidic glycans were chromatographically separated to obtain isomer-specific information and each glycan was relatively quantified by PGC-MRM MS. Total eleven *N*-glycan isomers were monitored.

Results

Glycans released from human serum were initially profiled by nanoLC/MS and about 100 glycan compositions were identified. Of them, sialylated glycans which are highly present in immunoglobulin-related proteins and in human serum were further selected for quantitative analysis. Total eleven *N*-glycan isomers consisting of three major compositions were successfully separated on a PGC column and quantitated using MRM-MS. In particular, four glycans including mono/di-sialylated bi-antennary *N*-glycans (Hex5 HexNAc4 NeuAc1-2) were highly expressed in BD patients compared with control group. These glycan isomers showed high diagnostic efficacy having AUC of ROC curve over 0.98. As well as disease-specific glycans, we found the glycan expression level of female control groups was ten times higher than that of male.

Conclusions

While global profiling of serum *N*-glycome has already revealed potential biomarkers for various diseases, glycan quantitation is analytically challenging due to the complexity caused by glycan isomers with a wide dynamic range of intensity. In this study we developed new technology to quantify isomer-specific glycans using PGC based MRM-MS for Behcet disease monitoring.

Novel Aspect

The first study of isomer-specific glycan quantitation using MRM MS for Behcet Disease Monitoring

References

1. S. Hua, et al., Analytical Chemistry, 85, 4636-4643 (2013)
2. H. J. An, et al., Journal of Proteome Research, 5, 1626-1635 (2006)

ThP-12 / GLYCOSAMINOGLYCAN ANALYSIS OF PROSTATE CANCER TISSUE MICROARRAYS

Gábor Tóth^(1,2) - **Oliver Ozohanics**⁽¹⁾ - **András Ács**⁽¹⁾ - **Ágnes Révész**⁽¹⁾ - **Károly Vékey**⁽¹⁾ - **László Drahos**⁽¹⁾ - **Lilla Turiák**⁽¹⁾

MS Proteomics Group, Research Center for Natural Sciences, Hungarian Academy of Sciences, Budapest, Hungary (1) - Department of Chemical Technology and Biotechnology, Budapest University of Technology and Economics, Budapest (2)

Keywords: glycosaminoglycan, HILIC-WAX chromatography, self-packed nanoHPLC column, tissue microarray, prostate cancer

Introduction

Cancer research is one of the most studied areas of science and prostate cancer (PCa) is one of the most common types of cancer among men. Tissue samples, especially biopsies are often used in mass spectrometry based biomarker research and have a great potential in understanding biochemical mechanisms underlying diseases such as cancer. Tissue microarrays (TMA) consist of several biopsies (generally 1.5 mm in diameter) of different patients fixed on a microscope slide.

Methods

Glycosaminoglycan (GAG) chains were digested into smaller fragments (GAG disaccharides) on the tissue surface with bacterial lyase enzymes prior to the nanoLC-MS analysis. The development and optimization of a chromatographic and MS method for the investigation of the GAG disaccharides was performed involving the use of novel self-packed HILIC-WAX capillary columns. For the statistical evaluation of datasets T-test and principal component analysis were used.

Results

The aim of our work was to develop and apply advanced nanoLC-MS techniques to reliably identify GAGs from the surfaces of PCa TMAs.

As for the separation method, effects of the ionic strength, the pH and the eluent strength were investigated and an isocratic separation method was found to be suitable for separating the GAG disaccharides. The limit of detection was ca. 1 fmol for each disaccharide and the limit of quantitation was between 10-50 fmol.

The nanoLC system was on-line coupled to mass spectrometer operated in negative mode and several TMA cores of different pathological grades of PCa biopsies were analyzed.

Statistically significant changes in heparan-sulphate composition were observed. The ratio of the different GAG disaccharides was found to be a potentially useful indicator of prostate cancer progression. A good correlation occurs between the cancer grade and sulphation patterns of the heparan-sulphate chains; the ratio of doubly and triply sulfated disaccharides increases with cancer progression.

Conclusions

GAGs play critical roles in cancer progression and our new isocratic HILIC-WAX-nanoHPLC-MS method provides extremely high sensitivity and reproducibility thus enabling the determination of the sulfation pattern from single TMA cores (<10 μ g tissue sample). The ratio of the different GAG disaccharides correlates with PCa progression; the amount of doubly and triply sulfated disaccharides shows an increase with cancer grade.

Novel Aspect

Our novel HILIC-WAX nanoHPLC-MS method enables determination of the GAG content of single PCa TMAs and sulfation pattern of these TMA cores may be a useful indicator of PCa grade. #

Acknowledgement

Gábor Tóth appreciates the financial support of the Mass Spectrometry Division of the Italian Chemical Society providing the IMSC 2018 participation fellowship.

ThP-13 / THIN-LAYER CHROMATOGRAPHY COUPLED WITH NANO-MATRIX BASED MALDI-TOF MASS SPECTROMETRY FOR STRUCTURAL ELUCIDATION OF GLYCANS

Elias Gizaw - Mei-Chun Tseng - Yuju Chen

Academia Sinica, Institute of Chemistry, Taipei

Keywords: Carbohydrate, Thin-layer chromatography, nanoparticle, MALDI-MS

Introduction

Our recent study demonstrated that 2,5-dihydroxybenzoic acid conjugated Fe₃O₄ nano-matrix (DHB@MNPs) was powerful for simultaneous ionization and fragmentation of oligosaccharides by MALDI MS [1]. However, analysis of more complex mixtures of carbohydrate requires further separation. In this study, we aimed to directly couple thin layer chromatography (TLC) separation of carbohydrate with DHB@MNPs based MALDI MS detection as a novel method of analysis.

Methods

Our new approach incorporates several stages of development. First, saccharides mixture separation was carried out on C18 modified reverse phase (RP) TLC plate. Second, we used ionic liquids as nanoparticle dispersion media. Third, dispersed DHB@MNPs was then deposited on TLC surface by spin coating technique. Fourth, the separated carbohydrates spots was directly scanned by MALDI TOF MS. Finally, the method was applied to human milk oligosaccharide analysis.

Results

Compared to normal phase (NP) TLC, carbohydrate samples were found to be better separated on RP-TLC and yields higher fragment ion intensity which may be due to weak retention of polar carbohydrate on non-polar RP-TLC. Spot homogeneity of spin coated and ionic liquid dispersed DHB@MNPs were evaluated by acquiring spectrum from different distance from the origin. Coefficient of variation (CV) of molecular and fragment ions intensity of each separated spot fall in acceptable range (CV < 15%). Besides the molecular ions of D-Glucose and β-Lactose, several fragment ions of Lewis Y were able to be detected by direct MALDI MS scanning which enabled unambiguous structural elucidation with single MS. The isomeric pairs of Lewis B and Y were also discriminated by DHB@MNPs-induced unique fragment ions [Z1α (372.3) and Y1α (390.4) unique for Lewis B; 2,4A3/Y2α (391.5), 0,3A3α (421.1) and Z2α (534.5) unique for Lewis Y]. Using our optimized method more than 32 neutral human milk oligosaccharides were identified from a donor at after one month lactation time.

Conclusions

With our integrated platform we were able to demonstrate separation and detection of a model saccharide mixtures. DHB@MNPs matrix was dispersed well by ionic liquid and uniformly deposited on TLC surface by spin coater. Automatic scanning allowing more accurate, faster and reproducible data acquisition. Overall, our approach provides alternative efficient carbohydrate analysis without the need of the common prior chemical derivatization and tandem MS.

Novel Aspect

RP-TLC separation and spin coating of ionic liquid dispersed DHB@MNPs followed by MALDI MS scanning as a facile method for structural elucidation of carbohydrate mixture.

References

1. Obena, R. P.; Tseng, M. C.; Primadona, I.; Hsiao, J.; Li, I. C.; Capangpangan, R. Y.; Lu, H. F.; Li, W. S.; Chao, I.; Lin, C. C.; Chen, Y. J., *Chemical science*, 6 (8), 4790-4800 (2015).

ThP-14 / GLYCOLIPIDOMICS OF HUMAN BRAIN BY COMBINED MODERN MASS SPECTROMETRY METHODS: CHIP NANO ESI MS AND ESI IMS MSⁿ

Mirela Sarbu⁽¹⁾ - *Željka Vukelić*⁽²⁾ - *David Clemmer*⁽³⁾ - *Alina Diana Zamfir*⁽¹⁾

National Institute for Research and Development in Electrochemistry and Condensed Matter, Department of Condensed Matter, Timisoara, Romania (1) - University of Zagreb Medical School, Department of

Keywords: chip nanoESI MS; ion mobility separation mass spectrometry; human adult hippocampus; gangliosides; collision-induced dissociation.

Introduction

The central nervous system possesses the highest content of gangliosides (GGs) [1]. Considered as biomarkers of brain development, aging, and certain diseases, the gangliosides expressed in an adult hippocampus (AH) specimen, the first brain region affected in neurodegeneration [2], were here systematically investigated in a combined approach based on chip nanoelectrospray ionization (ESI) QTOF MS and ESI ion mobility separation (IMS) MS.

Methods

GG solution of 5 pmol/ μ L concentration in methanol was investigated in negative ion mode using an orthogonal hybrid QTOF MS coupled to a fully automated chip-based nanoESI robot and in parallel by IMS using a Synapt G2s. A 2 kV ESI voltage was used for both experiments, while the sample flow rates were 50 nL/min for QTOF and 2 μ L/min for IMS. For the IMS experiments, IMS wave velocity was set at 650 m/s, IMS wave height at 40 V. CID MSⁿ were performed at energies within 10-50 eV.

Results

The generated data clearly showed similarities of ganglioside pattern. However, significantly lower extent of in-source-induced fragmentation of GGs using chip-based ESI versus the conventional ESI-based ion source was observed. This attribute represented the major advantage considering the content in polysialylated species and complex structures modified by fucosylation or acetylation, respectively, which are prone to in-source detachment. On the other hand, ESI IMS MS separated the components based on their differential mobility through a buffer gas, isolating the GGs based on the carbohydrate chain length and the degree of sialylation. Altogether, the optimized conditions for both experiments led to the identification of over 100 distinct ganglioside structures in human adult hippocampus, some of them modified by fucosylation, acetylation and *N*-acetyl galactosamine attachment. Furthermore, by applying CID MSⁿ, novel GG species, not identified before in AH, were structurally investigated in details.

Conclusions

In this work, we report on a systematic and comparative assay based on two high-performance MS platforms for the determination of GG structural diversity in a native mixture purified from normal adult human hippocampus. The results demonstrate that the combination of chip-based nanoESI QTOF MS and ESI IMS MS and MSⁿ data is proficient to unequivocally detect and characterize complex glycolipids with potential biomarker role and discover new structures.

Novel Aspect

The two-step approach based on chip- nanoESI QTOF MS and ESI IMS MS provided novel insights into the pattern and structure of GGs in adult hippocampus.

References

1. Sarbu M., Vukelić Ž., Clemmer D. E., Zamfir A. D., *Biochimie*, 139, 81-94 (2017).
2. Vukelić Ž., Zarei M., Peter-Katalinić J., Zamfir A. D., *J. Chromatogr. A*, 1130, 238-245 (2006).

ThP-15 / TOP DOWN PROTEOMICS OF HUMAN ALPHA-ACID-GLYCOPROTEIN: FAST AND UNAMBIGUOUS GLYCAN ANALYSIS

Melissa Bärenfänger - Bernd Meyer

University of Hamburg, Department of Organic Chemistry, Hamburg, Germany

Keywords: Top-Down Proteomics, Glycan Analysis, negative ion mode, Alpha-Acid-Glycoprotein, ESI-MS

Introduction

Glycosylation is an important posttranslational modification as it modulates biological functions of proteins. An altered glycosylation is found in a variety of diseases such as cancer, neuronal diseases and inflammation. Therefore, monitoring changes in glycosylation leads to a better understanding on how glycosylation influences the development and progression of these diseases and can be applied for biomarker development.

Methods

We present a top down approach for qualitative and quantitative analysis of highly glycosylated human alpha-acid-glycoprotein (AGP). After purification of AGP from human plasma, subsequent ESI-q/TOF-MS measurement in negative ion mode revealed the broad mass distribution caused by glycosylation. To support data interpretation, AGP was also analyzed after desialylation. In-house developed bioinformatics tools were used for data interpretation.

Results

Analyzing AGP as native glycoprotein and after treatment with neuraminidase reduced sample preparation and led to a fast and robust method for glycan analysis. We were not only able to determine the proteins glycosylation pattern but also to identify various genetic variants for AGP-1 and AGP-2. The variants did not show significant differences in their glycosylation pattern. The fully sialylated glycoprotein showed an almost complete termination with sialic acids of all antennae, carrying on average 16 to 17 sialic acid residues. To demonstrate the possible use of this method for biomarker development, we analyzed one sample from a patient diagnosed with liver cirrhosis and one diagnosed with hepatocellular carcinoma. In both cases, a drastically increased fucosylation over all glycosylation sites was observed.

Conclusions

Our method presents a way of detecting different genetic variants and determine the glycosylation pattern specific for each variant at the same time. Bioinformatics tools analyzed the complex heterogeneity in glycosylation unambiguously. Changes in glycosylation as *e.g.* increased fucosylation due to cancer diseases can be detected and quantified easily.

Novel Aspect

Ionizing highly sialylated glycoproteins in negative ion mode enhances signal intensity and can be applied for quantification of glycans in use as clinical biomarker.

ThP-16 / DEVELOPMENT OF A NEW METHOD FOR THE IDENTIFICATION OF FOLLICLE STIMULATING HORMONE (FSH) GLYCOFORMS INVOLVED HUMAN FERTILITY

Chiara Melchiorre - Anna Illiano - Gabriella Pinto - Andrea Carpentieri - Angela Amoresano

Department of Chemical Science, University of Naples Federico II, Italy

Keywords: FSH (Follicle Stimulating Hormone), human fertility, glycosylation, mass spectrometry

Introduction

Follicle-stimulating hormone (FSH) plays a key role in the proliferation and development of ovarian granulosa cells in women and spermatogenesis in men [1-3]. The glycosylation of FSH is extremely important since it has a significant influence on its biological activity. The quantitative determination of FSH at the various stages of the menstrual cycle is fundamental, for this reason, it is important to develop new methodologies that can lead to the simple and sensitive determination of FSH.

Methods

During preliminary experiments, standard FSH was digested with trypsin and the resulting peptide mixtures directly analyzed by LC-MS/MS methodologies. Specific peptides were selected on the basis of their molecular ion and their most intense fragments and used as targets to set up the Multiple Reaction Monitoring (MRM) method. The set of peptides, defined by their specific MRM transitions, can only be associated to the FSH protein leading to unambiguous identification and quantitation of FSH even within a complex biological matrix like serum with unmatched specificity and sensitivity.

Results

Tandem mass spectrometry analysis allows us to characterize the glycosylation profile of several FSH molecules including the pituitary form in young and postmenopausal women, and some existing commercial forms. The definition of the glycoforms was carried out using cutting-edge techniques such as ETD fragmentation using a high resolution Orbitrap mass spectrometer. An immunoaffinity chromatography protocol was performed on a serum sample to isolate selectively FSH protein and to identify it. Finally, an innovative methodology for its identification and quantification was set up by MRM mass spectrometry. For the optimization of the MRM method, more unique peptides were selected and monitored in order to uniquely identify FSH in complex matrices. The building up of a calibration curve of standard peptides has allowed to quantify FSH in human serum. The research activities are performed in collaboration with the MERCK Group, leader in FSH production.

Conclusions

This study is addressed at defining the relationship between the glycosylation pattern of FSH and its biological activity and defining similarities and differences between the various forms and the standard molecule. Because of the extreme importance in determining the exact amount of FSH during the various phase of the menstrual cycle as well as in therapy, the development of a sensitive and robust method for the identification and the quantification of specific FSH glycoforms by advanced mass spectrometry methodologies was of great importance.

Novel Aspect

The quantitative determination of FSH is typically measured with the immunological assays [4] that suffer from numerous limitations correlated with the presence of false positives. These limitations can be overcome by innovative mass spectrometry techniques in Multiple Reaction Monitoring mode.

References

1. A. Amoresano, R. Siciliano, S. Orru A. Sirna, R. Napoleoni, V. Altarocca, E. De Luca and P. Pucci. Eur. J. Biochem. 242, 608-618 (1996).
2. Emperaire JC, Review of Physiology (2015).
3. Daya S. Treatments in Endocrinology 3, 161-171, (2004).
4. Amin, H. K., and W. M. Hunter. Journal of Endocrinology 48.3: 307-317(1970).

ThPS-S02 - Probing Protein Structure and Dynamics

Chairs: Lars Konermann, Giancarlo la Marca

ThP-17 / STRUCTURAL DIVERSITY OF NUCLEOSOMES CHARACTERIZED BY NANOESI-MS AND STRUCTURAL CALCULATION

Kazumi Saikusa⁽¹⁾ - **Akihisa Osakabe**⁽²⁾ - **Daiki Kato**⁽²⁾ - **Sotaro Fuchigami**⁽³⁾ - **Aritaka Nagadoi**⁽⁴⁾ - **Yoshifumi Nishimura**⁽⁴⁾ - **Hitoshi Kurumizaka**⁽⁵⁾ - **Satoko Akashi**⁽⁴⁾

Hiroshima University, Grad. Sch. of Science, Hiroshima (1) - Waseda University, Grad. Sch. of Adv. Science

and Engineering, Tokyo (2) - Kyoto University, Grad. Sch. of Science, Kyoto (3) - Yokohama City University, Grad. Sch. of Med. Life Science, Yokohama (4) - University of Tokyo, IQB, Tokyo (5)

Keywords: histone tail, ion mobility mass spectrometry, theoretical structural calculation

Introduction

Nucleosome core particle (NCP) is the basic structural unit of nucleosomes and composed of four histone proteins and ca. 150 bp DNA strand. The behavior of histone tails in nucleosomes with a variety of lengths of DNA was investigated using nanoESI-MS and nanoESI-IM-MS. The MS data were evaluated by comparing the structures generated by computational calculation in solution and gas phases.

Methods

Positive-mode nanoESI-MS and nanoESI-IM-MS data were obtained using a Waters Synapt G2 HDMS mass spectrometer. Nucleosomes were prepared with different lengths of DNA strands, 147, 250, or 342 bp DNA. Samples for MS experiments were prepared with ammonium acetate solutions. Possible structures of NCP were investigated by theoretical calculation.

Results

We found that the charge-state distributions of the nucleosomes were classified into two groups, narrow and wide, which was supported by IM-MS. By systematic analysis with ESI-MS in the positive-ion mode, it was suggested that the range of charge state of nucleosomes is related to the length of free DNA region. When the DNA strand has a region that is not involved in the NCP core structural formation, charge-state distributions of nucleosomes are narrow. In contrast, when the extra region of DNA does not exist, the charge-state distribution is widely spread. This suggests that histone tail regions contribute to generate structural diversity of nucleosomes in the absence of free DNA region. This was also supported by theoretical structural calculation of NCP.

Conclusions

When eukaryotic chromatin is viewed under an electron microscope, NCPs are linked via histone-free DNA strands. Although further study is required, our findings that the free DNA regions affects the histone tails' behavior may provide biological insight into the regulation mechanism of eukaryotic transcription initiation.

Novel Aspect

It was suggested that flexibility of histone tail regions is affected by the existence of free DNA regions. [1]

References

1. Saikusa k. et al., submitted.

ThP-18 / NATIVE STATE ORGANIZATION OF OUTER MEMBRANE PORINS UNRAVELED BY HDX-MS

Danilo Donnarumma - Claudio Maestri - Pietro Ivan Giammarinaro - Luigi Capriotti - Erika Bartolini - Daniele Veggi - Roberto Petracca - Maria Scarselli - Nathalie Norais

GSK Vaccines, Siena, Italy

Keywords: HDx-MS; membrane proteins; OmpF; outer-membrane vesicles; structural mass spectrometry

Introduction

Hydrogen-deuterium exchange (HDX) associated with mass spectrometry (MS) is emerging as a powerful tool to provide conformational information about membrane proteins [1]. Unfortunately, as for X-ray and NMR, HDx performed on reconstituted in vitro systems might not reflect the in vivo environment. We took advantage of the outer-membrane vesicles (OMVs) naturally released by E. coli to carry out conformational

and dynamic analysis of native OmpF through HDx-MS.

Methods

A new protocol compatible with HDx analysis that avoids hindrance from the lipid contents was setup. Briefly OMVs purified from *E. coli* Δ ompA were precipitated with TCA and washed with ice-cold acetone. In order to minimize the back exchange, all sample preparation steps have been shortened as much as possible. Each step required less than 5 s and was followed by centrifugation for 4 min at 4°C.

Results

The extent of deuterium incorporation was in good agreement with the X-ray diffraction data of OmpF [2] as the buried β -barrels incorporated a low amount of deuterium, whereas the internal loop L3 and the external loops incorporated a higher amount of deuterium. Moreover, the kinetics of incorporation clearly highlights that peptides segregate well in two distinct groups based exclusively on a trimeric organization of OmpF in the membrane: peptides presenting fast kinetics of labeling are facing the complex surrounding environment, whereas those presenting slow kinetics are located in the buried core of the trimer.

Conclusions

The work presented here clearly demonstrates that HDx-MS is a versatile technique that can provide information on solvent accessibility and spatial arrangement of major bacterial outer-membrane proteins. This study was possible by taking advantage of the OMVs and by developing a protocol that allows the removal of most of the lipid content while limiting the back exchange.

Novel Aspect

Gaining information on accessibility and dynamics of antigens in native conditions will be of pivotal importance for the future rational design of efficient vaccines.

References

1. Pirrone, G. F. ; Iacob, R. E. ; Engen, J. R. *Anal. Chem.* 2015, 87 (1), 99– 118
2. Kefala G. ; Ahn C. ; Krupa M. ; Esquivies L. ; Maslennikov I. ; Kwiatkowski W. ; Choe S. *Protein Sci.* 2010, 19(5), 1117-25

ThP-19 / EFFECTS OF OXIDATIVE STRESS ON THE CONFORMATIONAL STABILITY OF HEMOGLOBIN MEASURED BY VARIABLE TEMPERATURE NESI-IMS-MS

Daniel Woodall - Tarick El-Baba - Christopher Brown - Daniel Fuller - Shannon Raab - David Clemmer

Indiana University, Indiana University, Bloomington

Keywords: protein stability, oxidative stress, ion mobility, PTMs, native MS

Introduction

Structural stability plays an integral role in regulating protein function and activity. Stability and protein-protein binding interactions can be regulated by post-translational modifications that affect subunit binding affinities of protein complexes. Here, we use ion mobility spectrometry-mass spectrometry (IMS-MS) to explore how oxidative modifications affect subunit binding affinity and conformational stability of the hemoglobin tetramer complex.

Methods

IMS-MS experiments were performed using a Waters SYNAPT G2 mass spectrometer. Hemoglobin solution was dialyzed against 150 mM NH_4OAc (pH 7.4), and diluted to 15 μM [tetramer] for analysis. Protein solution temperature was controlled using a custom-fabricated heated nanoelectrospray device, in which the temperature of analyte solution and nESI emitter can be heated. Collision cross sections (CCS) were

estimated using the protocol described by Ruotolo et al. [1]

Results

Unfolding and dissociation of the native hemoglobin tetramer is induced by incrementally increasing the solution temperature in the nESI source. Shifts in the average charge state, and CCS of the $\alpha_2\beta_2$ tetramer arise with increasing solution temperature, leading to dissociation into $\alpha\beta$ dimers, as well as populations of both holo and apo monomers. Stabilities of each oligomeric state and conformer are determined based on their midpoint dissociation temperature (T_d). The $\alpha_2\beta_2$ holotetramer has a $T_d = 58^\circ\text{C}$, while the $\alpha\beta$ holodimer has a $T_d = 62^\circ\text{C}$, indicating the dimer is more stable at elevated temperatures relative to the tetramer. A subpopulation of $\alpha\beta$ dimers exists where one of the monomer units is in the apo state ($\alpha h\beta a$). A significant portion of this subpopulation shows evidence of an oxidation PTM at βCys93 . βCys93 has been identified as a "hotspot" for oxidation, as well as having implications on β -subunit heme loss due to its proximity to the heme binding pocket. [2] The oxidized species has a midpoint $T_d = 57^\circ\text{C}$, while the unmodified species has a $T_d = 61^\circ\text{C}$.

Conclusions

The T_d of the oxidized $\alpha h\beta a$ dimer species is 5°C lower than the holodimer, indicating that βCys93 oxidation has a considerable destabilizing effect on the α/β subunit interface, even more so than the loss of the heme group. In addition to decreasing the subunit binding affinity, the oxidized species favor structures with more compact CCS at increased temperature, suggesting a structural collapse or compaction may be involved in destabilizing the dimeric units.

Novel Aspect

Thermal stability analysis by IMS-MS reveals conformational destabilization induced by oxidative stress in hemoglobin tetramer complex.

References

1. Ruotolo B. T., Benesch J. L., Sandercock A. M., Hyung S. J., Robinson C. V., Nature Protocols, 3, 1139 (2008).
2. Strader M. B., Kassa T., Meng F., Wood F. B., Hirsch R. H., Friedman J. M., Alayash A. I., Redox Biology, 8, 363 (2016).

ThPS-S03 - Natural substances and their complexes

Chairs: Pascal Gerbaux, Luisa Mattoli

ThP-20 / DEVELOPMENT OF LC-MS/MS ANALYTICAL PROTOCOL FOR ARTEMISININ AND ITS DERIVATIVE COMPOUNDS QUANTIFICATION IN THE ARTEMISIA ANNUA L. EXTRACTS FOR COSMETIC APPLICATIONS

Rosana M. Alberici⁽¹⁾ - ***Elaine C. Cabral***⁽¹⁾ - ***Gislaine C. Nogueira***⁽²⁾ - ***Renata Vardanega***⁽²⁾ - ***Maria Angela A. Meireles***⁽²⁾ - ***Marcos N. Eberlin***⁽¹⁾

Thomson Mass Spectrometry Laboratory, University of Campinas, Campinas (1) - LASEFI-School of Food Engineering, University of Campinas, Campinas (2)

Keywords: Artemisia annua, artemisinin, supercritical extraction, LC-MS/MS

Introduction

Artemisia annua L., a plant used as traditional Chinese herb medicines, is rich in sesquiterpene lactones

that present a wide range of biological activities, exhibiting therapeutic potential as anti-inflammatory and antimicrobial. Here, an improved liquid chromatography-tandem mass spectrometry (LC-MS/MS) protocol was developed to quantify artemisinin and its derivative compounds in *A. Annua* extracts obtained by supercritical CO₂ fractional extraction.

Methods

A. annua extracts, including essential oil, fraction rich in artemisinin and a residual fraction were obtained using supercritical CO₂ fractional extraction, under optimized temperature and pressure experimental conditions. Thin layer chromatography (TLC) and liquid chromatography-tandem mass spectrometry (LC-MS/MS) were used, respectively, for qualifying and quantifying, artemisinin and its derivative compounds in the extracts.

Results

The developed LC-MS/MS method covers both the unambiguous identification and simultaneous quantification of artemisinin and its derivative compounds, while TLC protocol was effective for artemisinin semi-quantification during the optimization of the supercritical extraction process, as a rapid screening method. Using multiple reaction monitoring (MRM) method, artemisinin, deoxyartemisinin and dihydroartemisinin were detected as ammoniated adducts. For validation purposes, an *A. annua* extract, denominated E4, obtained by supercritical extraction, was used. Calibration curves were prepared in the concentration range of 50-500 ng/mL. Analytical parameters such as linearity, precision, accuracy, limit of detection (LD) and limit of quantitation (LQ) were established. The data indicated that the linearities resulting from the three analytes of interest showed correlation coefficient ($r \geq 0.99$). LD and LQ were 25 and 50 ng/mL, respectively. The relative inter and intraday standard deviation (RSD) was <5%.

Conclusions

The developed analytical methodology can be used for quality control of *A. annua* extracts due to be able to quantify simultaneously artemisinin, deoxyartemisinin and dihydroartemisinin. The standardized extracts could be use in cosmetic formulations.

Novel Aspect

Standardized and certified extracts from *A. annua* using LC-MS/MS protocol could be used more safely in new and effective natural cosmetic products.

ThP-21 / METABOLOMIC UHPLC-QTOF METHOD TO ANALYZE TANNINS IN NATURAL COMPLEX PRODUCTS USING ALL-IONS FRAGMENTATION ACQUISITION

Michela Burico⁽¹⁾ - *Luisa Mattoli*⁽¹⁾ - *Pietro Traldi*⁽³⁾

Aboca Spa società agricola, Area Ricerca Fitochimica, Sansepolcro (1) - Città della Speranza, Istituto di Ricerca Pediatrica, Padova (2)

Keywords: Tannins, Mass Spectrometry, All-Ions, Metabolomic

Introduction

Tannins are a class of compounds naturally occurring in some plants. They can be classified into condensed tannins (or proanthocyanidins, PAs) and hydrolysable tannins (HTs), further sub-classified in galloyltannins (GTs) and ellagitannins (ETs). Structurally, they are the most complex class of polyphenols as they are polymeric molecules. Here is described a metabolomic LC-HRMS method using All-Ions Fragmentation acquisition (AIF)[1] [2] to increase natural tannins identification.

Methods

Green tea leaves dried-extract (GTE) has been analyzed by means of UHPLC-qToF using a reverse phase

chromatographic column, in high resolution (2 GHz) and using AIF. In negative ion mode, AIF has been acquired at four different energies: 0, 20, 30, 40 eV. The EICs of fragment ions, from the high energy scan data, have been correlated to the precursor ion, from the low energy scan data. The EICs of each tannin subclass typical-fragments has been used to get tannins semi-quantitative data.

Results

Tannins are polymers characterized by certain degree of polydispersity, so that actually it is not possible to find on the market reference standards for all the possible molecular structures. In this context, the analysis in AIF shows its potential to study natural complex products containing tannins, as using this algorithm the fragment ion of UHPLC eluted compounds is correlated to its precursor. Combining the MS/MS information with retention time, high resolution accurate mass and product/precursor ion intensity ratios, it has been possible the identification and annotation of various tannins in GTE. After the creation of a specific tannins in-house database and the study of the tannins fragmentation behavior, the following fragments common to each tannin subclasses have been identified: ions at m/z 125.0236 and m/z 169.0140 for GTs, m/z 300.9980 for ETs, m/z 289.0719 and m/z 577.1358 for PAs. The annotated GTE tannins have been semi-quantified by external regression curves of Gallic acid (ions at m/z 169.0140) for GTs, Ellagic acid (ions at m/z 300.9980) for ETs and Procyanidin B1 (ions at m/z 577.1358) for PAs.

Conclusions

If reference standards are available, tannin scan be analyzed by comparison with the corresponding compounds. In case this it is not possible, metabolomic screening by AIF is an effective alternative! The approach presented here has allowed us to identify more than twenty GTE tannins, increasing the quality of the information on tannins, opening new perspectives in routine quality control of tannins as well in understanding their fate in biological studies.

Novel Aspect

A new metabolomic method for analyzing tannins in natural complex products while Reference standards are unavailable on the market is now available.

References

1. Application note 5991-2465EN, Agilent Technologies Inc., (2013).
2. Naz S., Gallart-Ayala H., Reinke S. N., Mathon C., Blankley R., Chaleckis R., Wheelock C. E., Analytical Chemistry, 89, 7933-7942, (2017).

ThP-22 / CHARACTERIZATION OF BETA AMINO ACID ESTERS IN ONONIS SPECIES

Nóra Gampe - András Darcsi - Lilla Ladocsi - László Kursinszki - Szabolcs Béni

Semmelweis University, Department of Pharmacognosy, Budapest

Keywords: phytochemistry, beta amino acids, isoflavonoids, fragmentation, quantitative

Introduction

The non-protein amino acids are widespread in the flora as important molecules against insect herbivores. Based on our previous results, Ononis species contain beta amino acids mainly in the form of isoflavonoid glucoside esters. The beta amino acid structural motif is common in CNS drugs. As Ononis species are used in traditional medicine, we aimed to characterize these structures as well as to define their quantity.

Methods

For the structural identification HPLC-HR-MS/MS, authentic standards, isolation and NMR experiments were used. HPLC-MS was applied to screen for these molecules in various Leguminosae plants. Since the quantitation of the beta amino acid was performed for its free form, a prior hydrolysis method has been

developed. To minimize the matrix effect during the quantitative HPLC-MS experiments, a further SPE method was applied for purification purposes.

Results

Based on the MS/MS fragmentation of the investigated molecules, six isoflavonoid glucosides and two beta amino acids could be assumed, resulting in altogether 12 compounds. The structures of isoflavonoid moieties were identified as formononetin, pseudobaptigenin, medicarpin, maackiain, onogenin and sativanone because of their specific MS fragmentation pattern [1]. The HR-MS/MS results combined with NMR data and comparison with authentic standards proved that homopipecolic acid and homoproline esterify the isoflavonoid glucosides. The presence of the esters was investigated in various Leguminosae plants (*O. spinosa* and *arvensis*, *Trifolium repens*, *Glycine max*, *Medicago sativa*), however, those were only characteristic for *Ononis* species. The roots of *M. sativa* also contained trace quantities while none of the other plants contained these molecules at all. Homopipecolic acid content of *Ononis* species ranged from 1 mg to 35 mg/100 g dried plant, while homoproline levels were between 0.1 mg and 16 mg/100g.

Conclusions

The structures of 12 isoflavonoid beta amino acid esters were elucidated. Screening various Leguminosae plants, only *Ononis* species contained them in significant amount, thus, these molecules could serve as chemotaxonomic markers. Levels of the beta amino acids were also quantified in these plants. Comparing these values with the ones required to induce CNS effects [2], [3] one can conclude that the rational use of *Ononis* species as traditional remedy, is safe.

Novel Aspect

This work presents the structure of 12 new isoflavonoid beta amino acid esters and firstly describes the quantity of cyclic beta amino acids in *Ononis* species.

References

1. N. Gampe, A. Darcsi, S. Lohner, S. Béni, L. Kursinszki, *J. Pharm. Biomed. Anal.*, 123, 74–81 (2016).
2. A. Kragler, G. Höfner, K. T. Wanner, *Eur. J. Pharmacol.*, 519, 43–47 (2005).
3. T. Steffan, T. Renukappa-Gutke, G. Höfner, K. T. Wanner, *Bioorg. Med. Chem.*, 23, 1284–1306 (2015)

ThP-23 / INVESTIGATION OF PLANT EXTRACTS FOR MEDICAL RESEARCH

Manuel Gronbach⁽¹⁾ - ***Christina Oppermann***⁽¹⁾ - ***Jana Marielle Modenbach***⁽²⁾ - ***Matthias Sandler***⁽²⁾ - ***Markus Lerch***⁽²⁾ - ***Udo Kragl***⁽¹⁾

University of Rostock, Department of Chemistry, Rostock (1) - University Medicine Greifswald, Clinic and Polyclinic for Internal Medicine A, Greifswald (2)

Keywords: natural product extraction, HPLC, GC/MS, LC/MS, MSn

Introduction

Nowadays the study of chemical compounds from natural sources is an established method to obtain new active substances. This strategy is used for the treatment of many different diseases and natural products show anti-inflammatory, antidiabetic, analgesic and anticancer activity. [1] We engage ourselves with the extraction, separation and characterization of various natural materials to contribute to the treatment and prevention of intractable diseases.

Methods

We use a two step method for the extraction which includes the use of ultrasound and microwave radiation. The polar extracts are then separated via HPLC to obtain fractions containing multiple substances and pure compounds. The complex mixtures are analyzed with LC/MS and GC/MS. Full scan experiments

with two scan modes give us a good overview of the extract composition. For further information we apply target molecule scans and tandem mass spectrometry (MSn).

Results

Following our two step extraction method we were able to produce different plant extracts. We extracted various natural sources like sea buckthorn, chokeberry, cornelian cherry and Chinese medical plants using mainly ethanol and ultra pure water as extraction solvents. With the help of our LC/MS and GC/MS devices the identification of many compounds from different substance classes was possible. Especially vitamins, organic acids, isoflavones, anthocyanins and other flavonoids were often found. Different glycosilated forms of those structures were also present within the extracts. We managed to obtain different fractions of our plant extracts after method developing with our HPLC system. Both full extracts as well as defined substance mixtures have been applied in cell culture studies. Some of them showed promising results in reducing cellular injury. [2] Our ion trap MS device enabled us to perform some MSn experiments. That way we collected structural information on several unknown compounds within the extracts.

Conclusions

We were able to prepare plant extracts which then were separated into fractions containing defined substances mixtures. The cell biological studies showed promising results for some of the tested compounds and we will continue to further separate the mixtures. The question remains if the biological activity is a synergistic effect of multiple natural products or only one active substance. We also need to obtain the pure unknown compounds for a full characterization.

Novel Aspect

The natural materials we extract in the PePPP-project are applied in the study of hereditary liver and pancreas disease, namely Wilson's disease, pancreatitis and Johanson-Blizzard syndrome.

References

1. Chandrika P. U., Rao A. S., Sirisha Y., Chandra S. R., International Journal of Pharmacy and Pharmaceutical Research, 6, 399-410 (2016).
2. Oppermann C., Engel N., Ruth W., Nebe B., Kragl U., Journal of Chemical and Pharmaceutical Research, 7, 175-185 (2015).

ThP-24 / MS ANALYSIS OF PROTEIN INTERACTIONS RELEVANT FOR THE BIOSYNTHESIS OF APE

Kudratullah Karimi

Goethe University Frankfurt am Main, Institute of physical and theoretical chemistry, Frankfurt am Main

Keywords: Native-MS, natural products, protein interactions, multienzyme complex

Introduction

Polyketides are a large class of biologically active and structurally diverse compounds. They are synthesized by polyketide synthases (PKSs), classified as type I, II and III enzymes. The Aryl-polyene ester (APE) system is a PKS type II-like system, where ten unbound proteins are involved, but the exact process of the biosynthesis of APEs is mostly unknown [1,2]. Here we report our current results of protein interactions in the synthesis of APE.

Methods

To reconstitute the APE biosynthesis, all proteins that are involved in this process were purified from *X. doucetiae*. Components were mixed in vitro, and samples were desalted with Micro Spin Columns (Thermo scientific) into ammonium acetate or Tris buffer directly before MS analysis. ESI mass spectra of the enzymes were obtained in positive ion mode using a Synapt G2S (Waters) mass spectrometer.

Results

The acyl carrier protein (ACP) is a small protein in PKS systems, which is responsible for carrying substrates cargo between active sites [3]. We see binding of ACP to different proteins showing protein-specific differences in binding strength. Acyl-acyl carrier protein synthetase (AasS) is a monomer and forms dimers in presence of ACP, it binds one ACP in monomeric and two ACPs in dimeric form. Ketosynthase 2 (KS2) can oligomerize up to a tetramer, but the dominant species are the monomer and dimer; it binds ACP in the monomeric and dimeric states. We can show the regulating effect of ACP on the oligomerization states of the proteins. The thioesterase (TE) forms tetramers and is able to bind up to four ACPs. The ketoreductase (KR) forms tetramers, where the dimer is the building block for the tetramer. Surprisingly we find that in this system the ketosynthase 1 (KS1) and the chain length factor (CLF) as well as the two dehydratases (DH1/DH2) build heterodimeric complexes, and only in their heterodimeric form, they can interact with ACP.

Conclusions

Incubation of all the proteins of the PKS and MS analysis supports the hypothesis that the APE system is a PKS II-like system. ACP has a regulating effect on the oligomerization of the proteins. We could determine the oligomeric states and binding behavior of the proteins that are involved in the APE process. These results are the first step in reconstructing the process of the APE biosynthesis.

Novel Aspect

Native MS analysis of protein interactions in PKS II system with separate proteins that act together in a multienzyme complex, which are used iteratively to form polyketides.

References

1. T. A. Schöner, S. Gassel, A. Osawa, N. J. Tobias, Y. Okuno, Y. Sakakibara, K. Shindo, G. Sandmann, H. B. Bode, *ChemBioChem*, 17, 247–253, (2016).
2. P. Cimermancic, M. H. Medema, J. Claesen, K. Kurita, L. C. Wieland Brown, K. Mavrommatis, A. Pati, P. A. Godfrey, M. Koehrsen, J. Clardy, *Cell*, 158, 412–421, (2014).
3. J. Beld, K. Finzel, M. D. Burkart *Chemistry & Biology*, 21, 1293–1299, (2014).

ThP-25 / ANALYTICAL AUTHENTICATION OF OPHIOGOPON JAPONICUS SAMPLES: AN ESSENTIAL STEP TO SECURE ITS SUPPLY AND ITS USE FOR THE DEVELOPMENT OF A SAFE AND EFFECTIVE DERMO-COSMETIC INGREDIENT.

Stéphanie Richer⁽¹⁾ - **Delphine Laporte**⁽¹⁾ - **Laurie Verzeaux**⁽¹⁾ - **David Boudier**⁽¹⁾ - **Isabelle Fourasté**⁽²⁾ - **M-J Howes**⁽³⁾ - **MSJ Simmonds**⁽³⁾ - **Brigitte Closs**⁽¹⁾

SILAB, R&D Department, Brive la Gaillarde (1) - Laboratoire Pharmacophores Redox, Phytochimie et Radiobiologie, Toulouse (2) - Kew Garden, Royal Botanic Garden, Richmond (3)

Keywords: Ophiopogon japonicus, LC-MS/MS, Authentication, secure supply

Introduction

Ophiopogon japonicus tubers are used in Traditional Chinese Medicine and are currently traded for the treatment of pathological skins. Their supply has to be traced and secured to ensure safety and efficacy of the resulting dermo-cosmetic ingredient. They share botanical and morphological characteristics with other plant species, particularly Liriope spicata. To avoid falsifications, analytical authentication of these plants tubers is thus an essential step.

Methods

In this context, we performed a comparative analytical study of various samples of Ophiopogon japonicus tubers, from different supply regions. Phytochemical profiles obtained by ultra performance liquid

chromatography coupled with tandem high resolution mass spectrometry (UPLC-MS/MS) were then directly compared with those of *Ophiopogon japonicus* and *Liriopes picata* specimens from referent collections.

Results

Chromatograms obtained from referent samples revealed specific markers of *Ophiopogon japonicus* and *Liriopes picata*, in particular the range of homoisoflavonoids [1-5] and alkaloids. Identification of the species could be performed thanks to comparison of parent and fragment ions to components described in the literature.

Based on these data, our results demonstrated that one of the samples contains alkaloids not present in *Asparagaceae* (*i.e.* the family of *Ophiopogon japonicus* and *Liriopes picata*), thus suggesting falsification of this sample. Two other samples display trace amounts of homoisoflavonoids, a variation probably due to the cultivation region which does not totally refute *Ophiopogon japonicus*. The other samples have the specific markers of authentic specimen of *Ophiopogon japonicus*.

Conclusions

This analytical analysis of samples from different regions allowed us to authenticate samples of *Ophiopogon japonicus* comparatively to certified specimens. This approach will allow to secure future supplies of this starting raw material and to detect early some falsifications. This work is essential to guarantee safety and efficacy of the resulting dermo-cosmetic ingredient.

Novel Aspect

This study shows that an analytical approach is complementary to botanics to authenticate a plant species, avoid falsifications and secure its supply.

References

1. Wu Y., Dong Z., Wu H., Ding W., Zhao M., Shi Q., Wang Q., Food Research International, 57, 15-25 (2014).
2. Ye M., Guo D., Ye G., Huang G., American Society for Mass Spectrometry, 16, 234-243 (2005).
3. Lu X., Tong W., Wang S., Li J., Zheng J., Fan X., Liu L., Journal of Pharmaceutical and Biomedical Analysis, 134, 60-70 (2017).
4. Chen M. -H., Chen X. -J., Wang M., Lin L. -G., Wang Y. -T., Journal of Ethnopharmacology, 181, 193-213 (2016).
5. Shang Z. -P., Wang F., Zhang J. -Y., Wang Z. -J., Lu J. -Q., Wang H. -Y., Li N., Chinese Journal of Natural Medicines, 15, 11, 801-815 (2017).

ThP-26 / LC-ESI-MS/MS ANALYSIS OF SHORT-CHAIN PEPTIDES WITH ANTI-AGING PROPERTIES

Maftuna Shamshiddinova - Yong-Moon Lee

Chungbuk National University, Chungbuk National University/ College of Pharmacy, Cheongju, Republic of Korea

Keywords: short chain peptides, anti-aging, cosmetic formulations, LC-ESI-MS/MS analysis

Introduction

Skin aging is a natural event characterized by progressive structural and physiological changes in the skin. The main outcomes are the atrophy of the skin followed by loss of elasticity and a slowed metabolic activity [1]. Peptides like GHK, pal-KTTKS, and pal-GHK stimulate collagen renovation and shows anti-aging activity [2,3,4,5]. Thus, cosmetic producers interested in the introduction of peptides in cosmetic formulations is greatly increased in last years.

Methods

4 modified peptides: 2,5-dihydroxybenzoic acid-PLG-OH (Db-PLG), 2,5-dihydroxybenzoic acid-LG-OH (Db-LG), 2,5-dihydroxybenzoic acid-PG-OH (Db-PG) and 2,5-dihydroxybenzoic acid-GL-OH (Db-GL) was quantified in cosmetic samples using LC-ESI-MS/MS method with positive ion mode. The method was validated using quality parameters like linearity, calibration, correlation of determination, detection limit, precision and, recovery.

Results

Retention times for peptides were following: 18.29 min for Db-PLG, 22.3 min for Db-LG, 22.56 min for Db-GL and 6.74 for Db-PG. The calibration curves were found to be linear with an average correlation of determination R^2 to 0.9926 for Db-PLG, 0.9903 for Db-LG, 0.9986 for Db-PG and 0.996 for Db-GL.

Detection of peptides was performed with an AB SCIEX 3200 triple quadrupole system in the MRM and positive electrospray modes. MRM transitions for Db-PLG, Db-LG, Db-PG, and Db-GL were m/z 422/234, 325.1/250, 309.1/234 and 325.1/194.0, respectively.

Limits of quantification was established for a signal to noise ratio equal to 10. This corresponds to 40 ng for 2,5-dihydroxybenzoic acid-PLG-OH, 10 ng for 2,5-dihydroxybenzoic acid-LG-OH, 50 ng for 2,5-dihydroxybenzoic acid-PG-OH and 10 ng for 2,5-dihydroxybenzoic acid-GL-OH.

Recoveries were checked at three different concentrations and it was varied from 99.3% to 101.8%.

The method was applied to detect and quantitate tri- and di-peptides in different cosmetic formulations like skin, lotion, and gel.

Conclusions

LC-ESI-MS/MS method was developed for the quantification of tri- and dipeptides with anti-aging properties: Db-PLG, Db-LG, Db-PG, and Db-GL. The method was validated with the International Conference on Harmonization of Technical Requirements for Registration of Pharmaceutical for Human Use guidance note represents validation procedure of analytical methods. The method was applied to determine short chain peptides in cosmetic formulations.

Novel Aspect

newly synthesized short chain peptides in different cosmetic formulas were quantified by using LC-ESI-MS/MS method.

References

1. Delalle-Lozica N., Acta Clin Croat, 49, 529-536 (2010).
2. Maquart F. X., Pickart L., Laurent M., Gillery P., Monboisse J. -C., Borel J. -P., FEBS Letters, 238, 343-346 (1988)
3. Lupo M. P., Cole A. L., Dermatol and Therapy, 20, 343-349 (2007).
4. Lintner K., Peschard O., International Journal of Cosmetic Science, 22, 207-218 (2000).
5. Robinson L. R., Fitzgerald N. C., Doughty D. G., Dawes N. C., Berge C. A., Bisset D. L., International Journal of Cosmetic Science, 27, 155-160 (2005).

ThP-27 / FROM MASSIVE CHICKEN DEATH TO INSECT POISON DIET THROUGH THE MASS SPECTROMETRY PRISM

Elena Stashenko - Jairo Rene Martinez - Fausto Prada - Julian Lozano - Yesica Bautista - Jesus Zambrano - Nathalie Nãñez

Universidad Industrial de Santander, School of Chemistry, Bucaramanga, Colombia

Keywords: pyrrolizidine alkaloid crotalaria N-oxide monocrotaline

Introduction

Many chickens died because their food used sorghum contaminated with *Crotalaria retusa* (Fabaceae) seeds with high concentration of pyrrolizidine alkaloids (PAs), highly toxic substances. We determined the PAs and their N-oxides (PANOs), in *Crotalaria* genus plants growing in Colombia. Alkaloid distribution was determined in root, stem, leaves, flowers, and seeds of six *Crotalaria* species and in *Utethesia ornatrix* eggs, larvae and moths, in relation to herbivory.

Methods

U. ornatrix larvae consumed leaves of *C. retusa*, *C. pallida*, *C. spectabilis*, *C. incana*, *C. maypurensis*, and *C. nitens* maintained in an experimental garden. Solvent extraction, matrix solid-phase dispersion were used with plant material and insect parts. Pyrrolizidine alkaloids (PA) were determined both free and as their N-oxides. GC-MS analysis of PANOs required their previous reduction with Zn. Extracts were analyzed on an Exactive Plus orbitrap LC-MS system.

Results

High-resolution mass spectrometry afforded the exact molecular masses of both pyrrolizidine alkaloids and their N-oxides (PANO). Monocrotaline, creatonotine B, crispatine, integerrimine and their N-oxides were found both in plant material (0.01 – 21 mg/g) and insect specimens (0.01 – 1 mg/g). *U. ornatrix* eggs contained callimorphine N-oxide additionally. PA + PANO content in *Crotalaria* spp. plant parts varied according to the species, but was highest in seeds, followed by flowers, leaves, and roots. *C. retusa* showed the largest PA+PANO content among the species studied. *U. ornatrix* larvae preferred *Crotalaria* spp. unripe seeds (2 mg PA+PANO/g) and apical leaves (0.3 mg PA+PANO/g). They accumulate PA and PANO, which are transferred to moth and eggs. Alkaloid content was higher in eggs (6 mg/g) than in other *U. ornatrix* stages. PA/PANO ratios of 1200 were determined in larvae fed with *C. retusa* leaves. PA and PANO content in various *Crotalaria* spp. parts did not change significantly in response to herbivory by *U. ornatrix* or mechanical cutting.

Conclusions

PAs and PANOs were simultaneously determined with LC-MS. MSPD permitted alkaloid isolation from small sample amounts. *U. ornatrix* and *Crotalaria* spp. exemplify co-evolutionary adaptation. Pyrrolizidine alkaloids produced by the plant as a defense measure are used by *U. ornatrix* as a protection from predators. *Crotalaria* spp. seeds (29 mg/g) and *U. ornatrix* eggs (5 mg/g) showed the largest alkaloid content, in which the N-oxide form was prevalent (PANO/PA ~ 1200).

Novel Aspect

With good reproducibility, matrix solid-phase dispersion afforded plant and insect extracts with high alkaloid recovery.

ThP-28 / CHARACTERIZATION OF ENVIRONMENTAL MICROORGANISMS BY PROTEINS AND LIPIDS MALDI-TOF FINGERPRINTS

Vincent Guérineau - Morgane Barthélemie - Véronique Eparvier - David Touboul

CNRS-ICSN, Université Paris Saclay, Gif sur Yvette

Keywords: MALDI-TOF, microorganism, protein, lipid, database

Introduction

MALDI-TOF MS is an emerging technique for microbial identification, characterization and typing that has been successfully used in the clinical field. Though, studies related to the characterization of remain scarce. In order to fill this gap, protein profiles of various environmental microorganisms have been acquired and

classified. Moreover lipid profiles have been obtained and are expected to be complementary to protein profiles for rapid classification.

Methods

The lipids of each strain are extracted after maceration in an MTBE / MeOH mixture and then recovery of the organic phase and evaporation of the solvent. These lipid extracts are solubilized in a CHCl₃/MeOH mixture (2/1) and then analyzed by MALDI TOF in positive mode (DHB matrix at 20 mg/mL in THF) and negative (9-AA matrix at 10 mg/mL in MeOH). Protein extraction: (CH₃CN/H₂O/HCOOH: 50/15/35) then MALDI-TOF analysis (α -CHCA matrix at 10 mg/mL).

Results

We were able to confirm that the spectral fingerprint of a protein extract is dependent on the species of environmental microorganism, whether for a bacterium or a fungus. For example, 16S sequencing of strains belonging to genera such as Bacillus is not sufficiently discriminating for identification at the species level. On the other hand, clusters within the same genus appear well defined when observing the spectral imprints of these same strains. Unlike protein extracts, lipid fingerprints do not seem to depend on the species but only on the genus in most cases. Only data on fungi of the genus Fusarium allow discrimination at the species level. The study related to lipid fingerprints has shown that strain classification can be performed even though genetic sequencing has not been performed or if the protein fingerprint is too specific to be related to a known species.

Conclusions

In our study we observed a perfect concordance of "clustering" between genetics and spectral fingerprints. We have therefore been able to provide a proof of concept concerning the possible identification of the different environmental microorganisms by comparison of protein and / or lipid fingerprints and to show the important potential of MALDI-TOF mass spectrometry for the identification of these environmental microorganisms.

Novel Aspect

New MALDI-TOF methodology for the identification of environmental strains (phytopathogens, mutualist microorganisms, endophytes, etc.)

ThP-29 / MASS SPECTROMETRY TECHNIQUE FOR STUDYING THE PEPTIDES CYCLISATION BY NATIVE CHEMICAL LIGATION

Magdalena Wierzbicka⁽¹⁾ – *Olga Skwara*⁽¹⁾ - *Seçil Deniz Ozan*⁽²⁾ - *Piotr Stefanowicz*⁽¹⁾

Faculty of Chemistry, University of Wrocław, Wrocław, Poland (1) - Faculty of Science, Izmir Institute of Technology, Izmir, Turkey (2)

Keywords: Cyclopeptides, Native Chemical Ligation, SPPS, LC-MS, HR-MS

Introduction

Cyclopeptides, widely spread natural products, exhibit various biological activities. The recent reports indicate that the native chemical ligation allows the synthesis of cyclic peptides[1]. However, this process has complex mechanism and optimization of this synthetic method needs extensive analytical characterization of reacting system including mass spectrometric methods.

Methods

The precursors of cyclopeptides susceptible to transamidation reaction were synthesized according to the solid phase peptide synthesis procedure. The synthesis and cyclization progress were monitored in

microscale by Triple Quadruple LC-MS, and FT-ICR MS and MS/MS. The structure of protected cysteamine derivatives applied as substrates in this synthesis were confirmed by NMR spectroscopy, FT-ICR MS and MS/MS all with ESI ion source.

Results

The new cysteamine derivative protected with Mmt was synthesized according to the procedure developed by Barlos et al. [2], purified by LC and analyzed by ESI-FT-ICR-MS with MS/MS and NMR. During MS experiment the fragmentation in source was observed and highly stable Mmt⁺ ion appeared. Model peptides as well as natural ones were synthesized on solid phase. Then the N,S-transfer of acyl group, followed by transtioesterification was performed, which allowed the transformation of peptide into the form of MESNa thioester. The interaction of thioester moiety with N-terminal Cys results in cyclization of a peptide. The intermediates of this reactions were monitored by Triple Quadruple HPLC-MS as well as FT-ICR-MS and MS/MS. The application of proposed method for the synthesis of biologically active peptides will be discussed in presented poster.

Conclusions

The presented LC-MS techniques allows the monitoring of peptides metathesis including the formation of intermediates. The developed method will allow highly efficient cyclopeptides synthesis combined with their liberation from the solid support during SPPS. The combination of the native chemical ligation with N,S-acyl shift and transtioesterification will allow the synthesis of any desired cyclopeptide containing cysteine residue.

Novel Aspect

The novel aspect is an application of the mass spectrometric methods for real time monitoring and optimization of the complex, biologically important chemical reactions.

References

1. Taichi M., Hemu X., Qiu Y., P. Tam J. P., Org. Lett. 15, 2620-2623 (2013)
2. Barlos K., Gatos D., Hatzi O., N. Koch, S. Koutsogianni, Int. J. Peptide Protein Res. 47, 148-153, (1996)

ThP-30 / AN UNUSUAL TMS ARTIFACT IN EI MASS SPECTROMETRY - PERRHENIC ACID TMS

W. Gary Mallard - Yufang Zheng - Sanford Markey - James L. Little

NIST, MSDC, Gaithersburg, United States

Keywords: TMS derivatization artifacts, surface reaction

Introduction

Trimethyl silyl derivatives (TMS) are one of the most common methods for extending the range of compounds for which GC/MS methods can be used, however as with most methods where reactions are involved there are often undesirable by-products [1]. Here we report an unusual artifact that appears to be caused by the mass spectrometer filament – specifically Perrhenic acid-TMS (Re(O₃)-O-TMS) appears to be formed during the ionization process.

Methods

Pure samples of compounds were analyzed both as the native and TMS derivatized compounds. The derivatization was carried out with BSTFA and the analysis was done with an Agilent 5977 MS with an Agilent 8939 GC. Chromatography used 15M Agilent VF-5MS column, 40:1 split ratio, 270°C inlet temperature, program 50°C start, ramp 15°C/min to 300, ion source temperature 250K. All data were analyzed using AMDIS.

Results

Several TMS spectra had low intensity ions at m/z 307,309. In some spectra these ions were above the expected MW and in some they were unreasonable losses from the molecular ion. The AMDIS deconvolution showed that these ions did not have the same time history as the main ions in the spectrum – they always appeared slightly delayed. However, the ions could not be arising from a distinct compound formed in the derivatization process since the ion maxima were always nearly coincident with the derivatized compound regardless of the elution time. Some the data files had separations in time between the 307/309 maxima and that of other ions was sufficient that AMDIS could extract a distinct spectrum. Comparing this spectrum with the NIST library gave a match factor of over 900 to Perrhenic acid-TMS ($\text{Re}(\text{O}_4)\text{Si}(\text{CH}_3)_3$). Careful examination of the smaller peaks in the spectra when the Perrhenic acid spectrum could not be deconvolved showed that the other expected peaks (294/292, 279/277, and 263/261) were seen with the isotope ratios expected for Re.

Conclusions

TMS derivatives reacted with the rhenium filament in the ion source to produce Perrhenic acid-TMS. The reaction product was typically slightly delayed in time from the derivative, but in some cases the temporal separation was not sufficient to allow for deconvolution of the pure derivative spectrum. Aging of the filament over a 2 month period removed the active sites on the filament and the perrhenic acid TMS was no longer observed.

Novel Aspect

Perrhenic acid-TMS from reaction of the rhenium filament and TMS derivatives is reported. The time history was used to determine the identity using AMDIS.

References

1. J. L. Little, J. Chromatogr. A, 844 (1999) 1-22, updated version at https://littlesandsailing.files.wordpress.com/2011/05/silyl_2014_text.pdf

ThP-31 / COMPARING THE PROFILES OF YELLOW-LEGGED HORNET VENOM (VESPA VELUTINA NIGRITHORAX) TO EVALUATE THE INFLUENCES OF SEASONAL VARIATIONS

Thao Nhi Le⁽¹⁾ - David Da Silva⁽¹⁾ - Cyril Colas⁽²⁾ - Eric Darrouzet⁽³⁾ - Patrick Baril⁽⁴⁾ - Lucie Petit Leseurre⁽⁵⁾
- Benoît Maunit⁽¹⁾

Institut de Chimie Organique et Analytique - ICOA UMR7311, Université d'Orléans, Orléans, France (1) - Centre de Biophysique Moléculaire (CBM), UPR 4301, Institut de Chimie Organique et Analytique (ICOA), Université d'Orléans-CNRS, Orléans, France (2) - Institut de Recherche sur la Biologie de l'Insecte (IRBI), UMR CNRS 7261, Université de Tours, Tours, France (3) - Centre de Biophysique Moléculaire (CBM), UPR 4301, CNRS, Orléans, France (4) - CHIMEX, Groupe L'Oréal, Le Thillay, France (5)

Introduction

Our study aims to exploit the Yellow-legged hornet venom (*Vespa velutina nigrithorax*) in the cosmetic. It is known that environmental factors could conduce morphological and behavior modification of insect to adapt to climatic changes. Our hypothesis is that seasonal variation can have a qualitative and quantitative impact on venom. In general, insect venoms contain a complex mixture of molecules, including proteins, peptides and low molecular weight compounds.

Methods

V. velutina venoms from winter or summer hornets were separated into two fractions using 10 kDa filter (Sartorius). Protein profiles were achieved from >10 kDa fraction by gel electrophoresis, their density was measured by TLC Visualizer (CAMAG). In parallel, <10 kDa fractions were analyzed by two techniques on

UHPLC-qToF (MS or MS/MS) and MALDI-TOF MS to obtain peptide profiles. Their differences of seasonal samples were accessed by statistical analysis (PCA, PLSDA).

Results

The protein profiles obtained in gel present a mass distribution from 10-95 kDa, in which the most abundant bands were observed from 20-50 kDa. Density profiles achieved from TLC Visualizer at white light, however, do not reveal the significant difference between winter and summer samples in this mass range.

Besides, a separation method by UHPLC-qToF (MS and MS/MS) was developed that allowed us to analyze peptides and low molecular weight compounds in *V. velutina* venom. The rich peptide profile was observed in each sample. In parallel, to develop an analytical method for investigating the peptide profile on MALDI-TOF, four matrices were screened with venom. The peptide profiles obtained from MALDI-TOF correspond to one from LC-MS and in addition, they provide more information in the high mass range. For global comparison, the statistical analysis results for each analytical method are in process.

Conclusions

There is no significant discrimination between winter and summer protein profiles on gel electrophoresis. The methodology for analyzing the profile variation of venom samples issued from different periods was established. The PCA and PLSDA approaches are essential for evaluating the difference and the results are in process.

Novel Aspect

The intraspecific comparison of venom compositions will be carried out to achieve a completed view of Yellow-legged hornet venom.

ThP-32 / GC-MS/MS DETERMINATION OF CIS AND TRANS- Δ 9-THC ISOMERS IN CANNABIS SATIVA L

Stefano Valera - Gianna Allegrone - Federica Pollastro

Universita Piemonte Orientale, Dipartimento di Scienze del Farmaco, NOVARA, Italy

Keywords: Cannabis sativa L., cis- Δ 9-THC, hemp, GC-MS/MS.

Introduction

Hemp (*Cannabis sativa* L.) is a complex plant able to produce many chemical entities. Its biomedical relevance has been focused on Δ 9-tetrahydrocannabinol (Δ 9-THC) (1). Anyway Δ 9-THC exists in two geometrical isomers: cis and trans. The trans- Δ 9-THC is the active compound responsible for the narcotic action of hemp. The existence and the distribution of the cis isomer in cannabis is still unclear (2).

Aim of this study was the determination of the cis and trans- Δ 9-THC in cultivar of Cannabis belonging to different chemotypes.

Methods

A gas chromatography-tandem mass spectrometry (GC-MS/MS) method was developed for the characterization and quantification of cis and trans isomers of Δ 9-THC. Calibration curves, with internal standard (IS), were prepared from standard solutions of both molecules. Linearity, LOD and LOQ were evaluated.

Aerial parts (leaves and inflorescences) of 31 hemp varieties were extracted in acetone. The dry extracts, added of IS, were analyzed by the developed method.

Results

The method developed resulted suitable for the determination of both Δ 9-THC isomers.

The calibration curves were linear in the range 5-200 µg/mL, LOD and LOQ were assessed respectively 0.3 and 1.0 µg/mL.

The study demonstrated the presence of cis- Δ^9 -THC in all the varieties analyzed in a range of concentration of 0.003 mg/g and 1.41 mg/g. The trans- Δ^9 -THC isomer, always present in higher concentration of the cis isomer, was assessed in the range of concentration of 0.004 mg/g and 15.50 mg/g.

Conclusions

The cis- Δ^9 -THC was detected in all samples analysed and its amount evaluated in comparison with the more known trans- Δ^9 -THC isomer. The cis/trans rate is always less than 1 and decrease with the increment of the total amount of Δ^9 -THC present in the plant.

Novel Aspect

For the first time the cis- Δ^9 -THC isomer was characterized and determined in different cultivar of Cannabis sativa L. .

References

1. Hanuš L. O., Meyer S. M., Muñoz E., Scafati O. T. Appendino G. Natural Products Report 33, 1357-1392, 2016.
2. Smith R. M., Kempfert K. D. Phytochemistry 16, 1088-1089, 1977

ThPS-S04 - Environmental MS

Chairs: Encarnation Moyano, Cecilia Bergamini

ThP-33 / STRATEGY TO IDENTIFY INDUSTRIAL PEAK EMISSIONS BASED ON LONG-TERM HRMS DATA TIME SERIES

Sabine Anliker⁽¹⁾ - Martin Loos⁽²⁾ - Matthias Ruff⁽¹⁾ - Rahel Comte⁽¹⁾ - Kathrin Fenner⁽¹⁾ - Heinz Singer⁽¹⁾

Eawag, Environmental Chemistry, Dübendorf (1) - looscomputing, R&D, Zürich (2)

Keywords: High resolution mass spectrometry, non-target screening, short-term peak emissions, industrial wastewater, time profiles

Introduction

Micropollutant emissions from industrial wastewater are poorly understood. Challenges arise from missing substance information due to confidentiality of production data. Moreover, short-term emissions from industrial batch processes are unlikely captured by conventional grab sampling. Hence, we developed a strategy for comprehensive identification of potential industrial peak emissions from treated wastewater to surface water.

Methods

Daily effluent samples (24-h composite) were collected for three months at two wastewater treatment plants (WWTPs) receiving different proportions of industrial wastewater. The samples were then analyzed by large volume direct injection LC-HRMS (orbitrap technology) in positive and negative ESI mode. Data processing followed a non-target approach using the freely available enviMass workflow [1] to extract time profiles of all detected signals from full scan spectra.

Results

After data cleaning and cross-sample componentization (isotopologue and adduct grouping), more than 20'000 and 60'000 time profiles of a priori unknown compounds were detected at the two WWTPs. More profiles were detected at the WWTP treating larger amounts of industrial wastewater. Based on observed time profile characteristics of compounds with known origin we developed a workflow to differentiate between continuous domestic emission and potential industrial peak emission profiles. Of the considered descriptors, profile intensity spread was the most conclusive, *i.e.* values for peak emissions were twice as high as for domestic emissions. The percentage of profiles finally classified as potential industrial peak emissions differed substantially between the two investigated WWTPs, being 3-fold higher for the plant receiving more industrial wastewater. For several of the observed peak emissions the underlying compound was identified and the emission was further assigned to a pharmaceutical production site in the catchment of the WWTP.

Conclusions

Our strategy of source attribution based on time profile characteristics proved to correctly identify industrial peak emissions and accurately predicted a higher share of peak emission profiles for the WWTP treating larger volumes of industrial wastewater. This approach is currently being applied to a national field study to estimate the contribution of industrial discharges on the overall micropollutant contamination of Swiss surface waters via WWTPs.

Novel Aspect

This study presents a non-target approach for systematic inspection of entire HRMS full scan time series datasets to identify potential industrial peak emissions.

References

1. Loss M., enviMass v3.4, <http://www.looscomputing.ch/eng/enviMass/overview.htm>

ThP-34 / TRANSFORMATION PRODUCT FORMATION IN DRINKING WATER TREATMENT

Andrea M. Brunner - Dennis Vughs - Wolter Siegers - Cheryl Bertelkamp - Annemieke Kolkman - Thomas Ter Laak

KWR Watercycle Research Institute, Chemical Water Quality and Health, Nieuwegein, Netherlands

Keywords: Non-target screening, transformation products, mass spectrometry, drinking water, water technology

Introduction

Transformation products (TPs) are formed in the water cycle through biological and technological processes. Despite their potentially altered toxicity compared to parent compounds, TPs formed by drinking water treatment are not routinely monitored and remain elusive. This is mainly due to the technical challenges in analyzing the often unknown, low concentration compounds. Their analysis requires non-target HRMS/MS methods and novel data analysis approaches.

Methods

Here, we performed lab scale experiments to monitor TP formation of the three organic micropollutants carbamazepine, clofibric acid and metolachlor during rapid sand filtration and ozonation, two readily applied biotic and abiotic drinking water treatments, respectively. TP identification from non-target data was facilitated through prediction of potential TPs based on literature and models, halogenated and/or isotopically labeled parent compounds.

Results

The experimental results showed that the degradation of parent compounds did not per se lead to mineralization of the compound, but rather to an abundance of TPs, in often low concentrations. Some of these TPs were bigger and less polar than their parent compounds, which was somewhat unexpected. The identification of peaks representing TPs was straightforward and semi-automatic with the developed workflow, based on statistical testing and peak area filters. The suspect screening based on TP suspect lists manually curated from literature mining and prediction tools was efficient for TPs in the lists. However, the majority of TPs identified did not match suspect list entries. Furthermore, the structural identification of these features, as well as of isobaric suspects remained labor and time intensive.

Conclusions

The majority of TPs remained structurally unidentified, and for the majority of identified TPs toxicological risk assessment was missing. Follow-up work should target and hopefully alleviate these issues. Finally, the developed workflow can be applied to pilot-scale experiments to allow TP monitoring in actual drinking water production.

Novel Aspect

develop and test an efficient workflow to monitor TP formation and identify drinking water treatment specific TPs on a lab-scale

ThP-35 / MOLECULAR LEVEL INVESTIGATION OF OIL CONTAMINATED SOILS GENERATED DURING THE GULF WAR WITH TWO-DIMENSIONAL GAS CHROMATOGRAPHY(GCXGC) AND ULTRA-HIGH RESOLUTION MASS SPECTROMETRY

Eunji Cho⁽¹⁾ - *Moon-Hee Park*⁽²⁾ - *Young Hwan Kim*⁽²⁾ - *Guyoung Kang*⁽³⁾ - *Yunju Cho*⁽⁴⁾ - *Sunghwan Kim*⁽¹⁾

Kyungpook National University, Department of Chemistry, Daegu (1) - Korea Basic Science Institute, Biomedical Omics Group, Cheongju (2) - Hankuk University of Foreign Studies, Department of Environmental Sciences, Cheoin gu (3) - Green-Nano Materials Research Center, Department of Chemistry, Daegu (4)

Keywords: Gulf War, Oil Spill, Contaminated Soil analysis, Two dimensional GC, FT-ICR MS

Introduction

The oil spills occurred at the end of the Gulf War is one of the serious environmental problem in history. Although 30 years have gone after the gulf war, the world is still concerned about the environmental impact of spilled oil. It is important to characterize and understand the chemical compositions of oil spills for effective cleaning of the oil spill area. In this study, extracts from contaminated soil in oil spill area were analyzed using 2DGC and FT-ICR MS.

Methods

Two soil samples were collected each between 15~20cm below the surface (top soil) and about 1m below the surface (bottom soil). The soil samples were extracted using sonication in hexane and acetone. Soil extracts were filtered and dried. The extract and oil samples were further diluted with solvent for FT-ICR MS. For 2DGC analysis, Rtx-5MS (30m x 0.25mm, 0.25 μ m) is used as the first non-polar column and Rxi-17sil-MS (1.3m x 0.15mm, 0.15 μ m) as the second polar column.

Results

Class distributions obtained by (+) APPI, ESI and LDI MS did not show significant difference between the soil extracts and non-degraded crude oils. On the other hand, (-) ESI showed a noticeable difference. Increased abundance of O2 class was observed from the top soil extract. Relative abundance of O2 class compounds with DBE 2 or 3 were most abundant in the top soil extracts. And the A/C ratio(acyclic/ Σ (2-4 ring cyclic) O2 species) of topsoil extract was the lowest. [1] In addition, when the soil samples were extracted with water,

a significant number of polar components were detected only in the topsoil extract. About 300 peaks were observed from 2DGC MS data obtained from soil extracts and oil. Volatile alkane peaks were disappeared and only aromatic compounds were observed from soil extracts.

Conclusions

Comprehensive analyses of soil extracts and oil were done by use of 2D GC and FT-ICR MS. The most of the light compounds in the soil were evaporated because of hot weather in Kuwait. And the top soil extract showed significant abundance of oxidized peaks. However, it was surprising to find that extensive oxidation has not occurred in the bottom soil. This study shows that chemical compositions are different in contaminated soils located at different depth from the surface.

Novel Aspect

This is the first study to analyze contaminated soils from spilled oil during the Gulf War using 2DGC and Ultra-High resolution MS.

References

1. S Kim, LA. Stanford, RP Rodgers, AG Marshall, CC Walters, K Qian, LM Wenger, P Mankiewicz, *Organic Geochemistry*, 36, 1117-1134 (2005).

ThP-36 / OCCURRENCE OF ORGANOPHOSPHATE FLAME RETARDANTS IN SWIMMING POOLS AND THEIR ENDOCRINE DISRUPTING POTENTIAL

Pei-Hsin Chou - Wen-Ting Shih

National Cheng Kung University, Department of Environmental Engineering, Tainan City, Taiwan

Keywords: Organophosphate flame retardants, Endocrine disrupting activity, In vitro bioassays

Introduction

The use of organophosphate flame retardants (OPFRs) has been increasing because of the restriction and phase out of brominated flame retardants. Environmental contamination of OPFRs is of growing concern since they have been reported to show various toxic effects, such as carcinogenicity, neurotoxicity, and endocrine disrupting activity [1-4]. The aim of this study was to investigate the occurrence of OPFRs in swimming pools and their potential health effects.

Methods

Swimming pool water samples were separated into dissolved phase (W) and suspended solids phase (SS) by filtration. W and SS samples were extracted by solid phase extraction and Soxhlet extraction, respectively. Concentrations of selected OPFRs were quantified using an isotope-dilution liquid chromatography-tandem mass spectrometry (LC-MS/MS) method. Endocrine disrupting activities of OPFRs and pool samples were investigated using in vitro bioassays [5].

Results

LC-MS/MS results revealed that 7 and 6 OPFRs selected for analysis were detected in W and SS samples, respectively, and the detection frequencies ranged from 65 to 100%. Tris(1-chloro-2-propyl) phosphate was the dominant OPFR in W samples, which was the major OPFR found in Australian public swimming pools as well [6]. Tri-n-butyl phosphate was the dominant OPFR in SS samples, and its concentrations were about 3 orders higher than other OPFRs. Bioassay results showed that selected OPFRs exhibited antagonist activities for androgen receptor (AR), glucocorticoid receptor (GR), and progesterone receptor (PR). AR and GR antagonist activities were also frequently found in W samples collected from southern Taiwan. By contrast, no AR antagonist activities were found in SS samples, and PR agonist/antagonist activities were not detected in W or SS samples as well.

Conclusions

The present study demonstrated that OPFRs were ubiquitously present in W and SS samples collected from Taiwanese swimming pools. Also, OPFRs exhibited AR, GR, and PR antagonist activities, suggesting potential health risks of OPFR exposure from swimming pool water. Further investigation will be undertaken to evaluate the contribution of OPFRs to the receptor antagonist activities detected in pool samples.

Novel Aspect

Trace concentrations of OPFRs in Taiwanese swimming pools were successfully analyzed for the first time by using the isotope-dilution LC-MS/MS method.

References

1. van der Veen I., de Boer J. V., Chemosphere, 88, 1119-1153(2012).
2. Chen G., Jin Y., Wu Y., Liu L., Fu Z., Environmental Toxicology and Pharmacology, 40, 310-318(2015).
3. Kojima H., Takeuchi S., Itoh T., Iida M., Kobayashi S., Yoshida T., Toxicology, 314, 76-83(2013).
4. Zhang Q., Wang J., Zhu J., Liu J., Zhao, M., Environmental Science & Technology, 51, 5803-5810 (2017).
5. Ito-Harashima S., Shiizaki K., Kawanishi M., Kakiuchi K., Onishi K., Yamaji R., Yagi T., Journal of Pharmacological and Toxicological Methods, 74, 41-52 (2015).
6. Teo T. L., Coleman H. M., Khan S. J., Science of the Total Environment, 569-570, 469-475 (2016).

ThP-37 / ASSESSMENT OF BPA AND DEHP METABOLITES IN URINE SAMPLES FROM MOTHER-CHILD PAIRS OF LIFE PERSUADED PROJECT

Veronica Della Latta⁽¹⁾ - **Fabrizia Carli**⁽¹⁾ - **Emma Buzzigoli**⁽¹⁾ - **Demetrio Ciociaro**⁽¹⁾ - **Melania Gaggini**⁽¹⁾ - **Anna Paola Pala**⁽¹⁾ - **Andrea Raffaelli**⁽¹⁾ - **Raffaele Conte**⁽¹⁾ - **Luca Busani**⁽²⁾ - **Francesca Maranghi**⁽²⁾ - **Sabrina Tait**⁽²⁾ - **Stefano Cianfarani**⁽³⁾ - **Cinzia La Rocca**⁽²⁾ - **Amalia Gastaldelli**⁽¹⁾

National Research Council, CNR, Institute of Clinical Physiology, IFC, Pisa (1) - Istituto Superiore di Sanità, Toxicology Unit, Roma (2) - Ospedale Pediatrico Bambino Gesù, Diabetologia, Roma (3)

Keywords: Ultra-High Pressure Liquid Chromatography Quadrupole Time of Flight (UHPLC/MS-QTOF), di-2-ethylhexyl phthalate (DEHP), bisphenol A (BPA), Gas Chromatography-Mass Spectrometry (GC-MS)

Introduction

BPA and DEHP are endocrine disruptors (ED) widely used as plasticizers and associated with the development of endocrine and metabolic diseases. In the LIFE PERSUADED project we developed UHPLC/MS-QTOF and GC-MS methods to monitor exposition to BPA and DEHP (by the evaluation of metabolites MEHP, MEHHP, MEOHP) in 900 healthy pairs mother-child living in rural and urban area from North to South Italy. Questionnaires on lifestyle and nutrition were recorded[1-3].

Methods

Urine samples(500 µl) were treated with enzymatic deconjugation to hydrolyze the glucuronide and separate free phthalates and BPA. Compounds with ¹³C standard were purified with C18SPE. DEHP metabolites were quantified by UHPLC/MS-QTOF (ZORBAX SB Phenyl 2.1x100mm 1.8 µm); the sample was then derivatized with BSTFA 1%TMS and acetonitrile and BPA was quantified by GC-MS in selective ion monitoring(J&W DB-5MS 30m, 0.25mm, 0.25µm) and normalized by creatinine.

Results

We were able to quantify DEHP metabolites in all urine samples. Preliminary data in 350 children (4-14 years, males, M and females, F) and mothers (Mo) showed that range were: MEHP 1.3-187.5, 1.49-117.26 and 1.9-51.8; MEHHP 0.74-1502.68, 0.0052-953.75 and 2.69-165.3; MEOHP 0.82-303.82, 1.27-827.03 and 0.024-51.96µg/g of creatinine (respectively in F, M vs Mo).

BPA analyses by UHPLC/MS-QTOF were affected by a matrix effect that was evident only in deconjugated samples. For this reason, after UHPLC/MS-QTOF analysis samples were dried and derivatized with BSTFA to be analyzed in GC-MS. Total BPA range was 0.69-1036.55 µg/g of creatinine(F) and 0.951-1146.786 µg/g of creatinine(M) in children and 0.84-43.78 µg/g of creatinine in their mothers.

Internal exposure was different according to living area, age (young children showed higher values), with gender differences. Data evaluation of questionnaires and food diaries is in progress and will be integrated with ED exposure data to provide the determinants to be used in risk assessment.

Conclusions

Although not persistent, BPA and DEHP metabolites have been detected in all enrolled children, with different exposure was different according to living area, age (lower range), with gender differences.

Novel Aspect

We set up a method to quantify DEHP metabolites by UHPLC/MS-QTOF and BPA concentration by GC-MS starting from 500 µl of urine.

References

1. Frederiksen H., Skakkebaek N. E., Andersson A.-M., Mol. Nutr. Food Res., 51, 899-911 (2007).
2. Covaci A., Den Hond E., Geens T., Govarts E., Koppen G., Frederiksen H., Knudsen L.E., Mørck T.A., Gutleb A.C., Guignard C., Cocco E., Horvat M., Heath E., Kosjek T., Mazej D., Tratnik J.S., Castaño A., Esteban M., Cutanda F., Ramos J.J., Berglund M., Larsson K., Jönsson B.A., Biot P., Casteleyn L., Joas R., Joas A., Bloemen L., Sepai O., Exley K., Schoeters G., Angerer J., Kolossa-Gehring M., Fiddicke U., Aerts D., Koch H.M. Environ Res., 141, 77-85 (2015).
3. Černá M., Malý M., Rudnai P., Középesy S., Náray M., Halzlová K., Jajcaj M., Grafnetterová A., Krsková A., Antošová D., Forysová K., Den Hond E., Schoeters G., Joas R., Casteleyn L., Joas A., Biot P., Aerts D., Angerer J., Bloemen L., Castaño A., Esteban M., Koch H.M., Kolossa-Gehring M., Gutleb A.C., Pavloušková J., Vrbík K. Environ Res., 141, 118-24 (2015).

ThP-38 / DULOXETINE: DEGRADATION PROCESSES, TRANSFORMATION PRODUCTS ELUCIDATION BY MASS SPECTROMETRY AND *IN SILICO* ANALYSIS

Rodrigo Osawa⁽¹⁾ - **Maria da Conceição Oliveira**⁽²⁾ - **Olinda Monteiro**⁽¹⁾ - **Maria Helena Florêncio**⁽¹⁾

Faculdade de Ciências, Universidade de Lisboa, Lisboa (1) - Instituto Superior Técnico, Universidade de Lisboa, Lisboa (2)

Keywords: Time-of-flight mass spectrometry; duloxetine; transformation products, *in silico* toxicity prediction

Introduction

This study aimed to elucidate the transformation products (TPs) of duloxetine (DUL), a serotonin-norepinephrine reuptake inhibitor, under controlled experimental conditions, to simulate processes that can occur in the environment and WWTP: photodegradation under UV irradiation and chlorination. Moreover, the ecotoxicity of the transformation products was predicted based on *in silico* tools [1]

Methods

The experiments were carried out on spiked distilled water. The chlorination consisted of adding sodium hypochlorite in the solution and for the photodegradation an UV irradiation source was used (450 W medium-pressure mercury-vapor lamp). The time profiles of formation/degradation of the transformation products were analyzed. Samples were directly injected in the LC-QTOF/MS system in full-scan and MS/MS mode.

Results

The parent compound completely degraded in the photodegradation experiments after 30 minutes of UV irradiation whereas the chlorination process accounted for about 90% of its removal. On the whole, eleven TP's were formed, nine from the contact with UV light and two from the chlorine reactions. These TP's were identified and their structures elucidated by analysis of their fragmentation patterns. Their formation has been attributed to, namely, oxidation, hydrolysis and dehydrogenation reactions. All TP's formed were totally degraded in both experiments: photodegradation in 4 hours and chlorination in 1 day. In comparison with DUL all TP's presented low ecotoxicity with the endpoints LC₅₀ F. minnow, *D. magna* and *T. pyriformis* from ECOSAR tool, and the quantitative structure–activity relationship (QSAR-TEST) analysis for mutagenicity resulted in five TP's with Ames-positive prediction.

Conclusions

Our experiments showed that duloxetine reacts with chlorine and under UV irradiation, producing transformation products, the structures of which were elucidated by analysis of their fragmentation patterns obtained using a LC-QTOF/MS equipment. These experiments also showed that the TP's and the parent compound degraded in both processes. Moreover the in silico analysis enabled to conclude that some TP's presented mutagenic potential.

Novel Aspect

Duloxetine is frequently detected in effluents, but little is known about its TP's. This study enables to throw some light on the behavior of this compound in the environment.

References

1. Martin, T. Toxicity Estimation Software Tool (TEST). U. S. Environmental Protection Agency, Washington, DC 2016

ThP-39 / MULTI-MODE IONIZATION TECHNIQUES FOR EFFICIENT MS/MS ANALYSIS OF PESTICIDES

Heather Gamble⁽¹⁾ - **Charles Jolliffe**⁽¹⁾ - **Avinash Dalmia**⁽²⁾

Perkin Elmer, R&D, Woodbridge (1) - Perkin Elmer, R&D, Shelton (2)

Keywords: MS/MS, multiple ionization modes

Introduction

Triple quadrupole mass spectrometry is well suited for targeted quantitative analysis of large numbers of compounds including pesticides, which can belong to a wide range of chemical classes and span a wide range of volatilities. Different pesticides tend to respond to a variety of ionization and separation techniques, such as LC-ESI, LC-APCI, GC-EI, etc. It is advantageous to support analysis of as many of these compounds as possible on one hardware platform.

Methods

A test suite of pesticides spanning different chemical classes and volatilities was selected for analysis using a single hardware platform consisting of a QSight 300 triple quadrupole mass spectrometer fitted with a versatile dual ion source capable of analyzing LC and GC amenable compounds using a variety of ambient ionization techniques.

Methods were developed to demonstrate highly sensitive analysis of diverse pesticides using a common hardware platform.

Results

Using a highly sensitive QSight 300 triple quadrupole mass spectrometer equipped with a versatile dual ion source it was possible to efficiently and quantitatively analyze both LC and GC amenable pesticides,

including difficult to ionize compounds such as chlordane or pentachloronitrobenzene. Multiple strategies are possible. One approach is to use HPLC to separate all compounds. The output of the LC was connected by T or switching valve to both probes on the dual ion source which was fitted with one LC-ESI and one LC-APCI probe. Traditionally difficult to ionize compounds, often analyzed by GC, can respond to APCI. Fast mode switching between the ESI and APCI probes was used to analyze both LC and GC amenable compounds using a single method, and within a single injection. The dual probe ion source can also be coupled to two separation devices, e.g. two LC's, or one LC, one GC, allowing for additional versatility.

Conclusions

A Q-Sight 300 triple quadrupole mass spectrometer equipped with a versatile dual probe source was used to demonstrate efficient, highly sensitive, quantitative analysis of both LC and GC amenable pesticides using a single hardware platform.

Novel Aspect

Simultaneous analysis of both LC and GC amenable pesticides using a single hardware platform.

ThP-40 / MULTICLASS ANALYSIS OF ANTIBIOTICS IN SEDIMENTS BY LC-Q-ORBITRAP HRMS. APPLICATION TO THE NERA RIVER

Danilo Giusepponi⁽¹⁾ - *Fausto Scoppetta*⁽²⁾ - *Marta Sargenti*⁽²⁾ - *Giorgio Saluti*⁽¹⁾ - *Simone Moretti*⁽¹⁾ - *Rosanna Rossi*⁽¹⁾ - *Marinella Capuccella*⁽²⁾ - *Roberta Galarini*⁽¹⁾

Istituto Zooprofilattico Sperimentale dell'Umbria e delle Marche "Togo Rosati", Centre for Development and Validation of Methods, Perugia (1) - Istituto Zooprofilattico Sperimentale dell'Umbria e delle Marche "Togo Rosati", Umbria Regional Center of Veterinary Pharmacovigilance, Perugia (2)

Keywords: veterinary drugs, sediments, liquid chromatography-high resolution mass spectrometry, Nera River

Introduction

Environmental research showed that pharmaceuticals from wastewater treatment plants, industries and hospitals can contaminate waters, sediments, aquatic flora and fauna [1]. Aquaculture plants also play an important role in this contamination. This study focused on the analysis of 58 antibiotics in sediments of Nera river (Umbria, Italy), an area in which several trout farms are present.

Methods

Sample preparation protocol consisted in three consecutive extraction steps followed by a tandem SPE clean-up (Oasis HLB + Strata-X-C). The LC-Q-Orbitrap conditions were the same described in our previous work [2]. The validation study has been performed at 1 and 10 µg/kg dry weight (dw). Samples from Nera River were collected monthly at three sites from March to December 2017.

Results

The optimized sample treatment involved three extractions with mixtures at different pHs (7, 3 and 10, respectively) to take out the wide range of substances included in the method scope: 58 analytes belonging to ten different antibiotic classes widely used in veterinary practices. After validation, the sediments collected from Nera River were analyzed (30 samples). Two tetracyclines (oxytetracycline and tetracycline), three quinolones (ciprofloxacin, flumequine and oxolinic acid) and trimethoprim were measured. In particular, oxytetracycline (and its epimer) was detected in the 100% of the samples with concentrations from 2 to 159 µg/kg dw. The other drugs were always lower than 10 µg/kg dw and their detection frequencies ranging from 3% (tetracycline) to 73 % (ciprofloxacin).

Conclusions

The results revealed the presence of drugs used in trout farming, but other antibiotics were from other sources (ciprofloxacin and oxolinic acid). Nera River collected at the same time [3]. On the other hand, it is well known that the sediment composition favors the adsorption of only certain drugs.

Novel Aspect

This is the first survey on the presence of antibiotics in the Nera River sediments.

References

1. Obimakinde S., Fatoki O., Opeolu B., Olatunji O. Environmental Science and Pollution Research, 24, 3274-3297 (2017).
2. Moretti S., Dusi G., Giusepponi D., Pellicciotti S., Rossi R., Saluti G., Cruciani G., Galarini R. Journal of Chromatography A, 1429, 175-188 (2016).
3. Giusepponi D., Moretti S., Saluti G., Scoppetta F., Capuccella M., Galarini R., Book of Abstract - 8th Symposium on Recent Advances in Food Analysis, p. 405 - 7-11 November 2017 - Prague (CZ).

ThP-41 / DETERMINATION OF 16 POLYCYCLIC AROMATIC HYDROCARBONS BY CORE-SHELL MOLECULARLY IMPRINTED POLYMERS PREPARED BY REVERSIBLE ADDITION FRAGMENTATION CHAIN TRANSFER POLYMERIZATION WITH APGC-MS/MS

Ali-Asghar Golbabanezhad-Azizi - Fereshteh Shahhoseini - Christina Bottaro

Memorial University of Newfoundland, Department of Chemistry, St. John's

Keywords: Magnetic molecularly imprinted polymers; dispersive solid phase extraction; atmospheric pressure chemical ionization gas chromatography; reversible addition fragmentation chain transfer polymerization

Introduction

Polycyclic aromatic hydrocarbons (PAHs) with carcinogenic and mutagenic properties are considered priority pollutants by regulatory agencies world-wide. Liquid-liquid extraction and solid phase extraction (SPE) are non-selective techniques influencing the sensitivity, accuracy and precision of MS analysis [1]. Molecularly imprinted polymers (MIPs) as selective sorbent can improve the efficiency, and reproducibility, and reduce matrix effects [2].

Methods

MIP particles were prepared by surface polymerization onto magnetic $\text{Fe}_3\text{O}_4@\text{SiO}_2$ nanoparticles through reversible addition fragmentation chain transfer (RAFT) [3]. The RAFT-MIPs were used for dispersive SPE (dSPE) of 16 PAHs from aqueous matrices, then isolated by magnetic collection. The analytes were desorbed into organic solvent prior to quantitation with atmospheric pressure chemical ionization gas chromatography tandem mass spectrometry (APGC-MS/MS).

Results

To obtain maximum performance by the RAFT-MIPs, the composition of polymeric coating as well as the parameters influencing extraction efficiency for enrichment of PAHs were assessed and optimized. The synthesized particles were also analyzed using XRD and IR measurements to confirm formation of selective binding sites through surface polymerization onto core-shell nanoparticles. Under optimized conditions (15 mg RAFT-MIPs, 40-mL sample volume, 1 min ultrasonic dispersion extraction time, 1 mL desorption solvent and 10 min desorption time), RAFT-MIPs were used for selective extraction of 16 PAHs from water samples. Good recognition properties (IF between 1.7 and 2.7) for extraction of PAHs in water samples were obtained without matrix manipulation, such as filtration or removal of interferents. The developed method with good detection limits (ng/L levels) and high precision (RSD<8%) provides a sensitive and reliable

procedure for analysis of PAHs in any complex matrix.

Conclusions

RAFT polymerization, as a controllable technique relying on living radical initiators, was used to create cavities with accessible and homogenous binding sites for the target molecules. The RAFT-MIPs yielded a rapid, selective and efficient procedure for extraction of PAHs from water samples. Indeed, application of a hydrophilic comonomers enhancing the water compatibility of magnetic MIP makes the sorbent suitable for the analysis in aqueous matrices.

Novel Aspect

Water compatible magnetic RAFT-MIPs were produced for rapid, selective and simultaneous extraction of 16 PAHs with reduced interferences for instrumental analysis.

References

1. Liew C. S. M., Li X., Lee H. K., Analytical Chemistry 88, 9095-9102 (2016).
2. Jiang T., Zhao L., Chu B., Feng Q., Yan W., Lin J. M., Talanta, 78,442-447 (2009).
3. Wu X., Wang X., Lu W., Wang X., Li J., You H., Chen, L., Journal of Chromatography A, 1435,30-38 (2016).

ThP-42 / FATE AND BEHAVIOR OF OIL SANDS NAPHTHENIC ACIDS IN A PILOT-SCALE TREATMENT WETLAND AS CHARACTERIZED BY NEGATIVE-ION ELECTROSPRAY IONIZATION ORBITRAP MASS SPECTROMETRY

Chukwuemeka Ajaero⁽¹⁾ - *Kerry Peru*⁽²⁾ - *Monique Simair*⁽³⁾ - *Vanessa Friesen*⁽³⁾ - *Sarah Hughes*⁽⁴⁾ - *Dena McMartin*⁽⁵⁾ - *Gwen O'Sullivan*⁽⁶⁾ - *John Headley*⁽²⁾

Environmental Systems Engineering, University of Regina, Regina, Saskatchewan, 3737 Wascana Parkway, S4S 0A2, Canada. (1) - Watershed Hydrology & Ecology Research Division, Water Science & Technology Directorate, Environment & Climate Change Canada, 11 Innovation Boulevard, Saskatoon, Saskatchewan, S7N 3H5, Canada (2) - Contango Strategies Limited 15-410 Downey Road, Saskatoon, Saskatchewan, S7N 4N1, Canada. (3) - Shell Health, Americas, One Shell Plaza, 910 Louisiana, Houston, TX, 77002, USA. (4) - Department of Civil, Geological and Environmental Engineering, 57 Campus Drive, University of Saskatchewan, Saskatoon, Saskatchewan, S7N 5A9, Canada (5) - Department of Earth & Environmental Science Mount Royal University 4825 Mount Royal Gate SW Calgary, AB, Canada, T3E 6K6 (6)

Keywords: Double bond equivalent, Naphthenic acids, Constructed wetland, Kinetic analysis, Orbitrap-MS

Introduction

Oil sand process-affected water (OSPW) is a complex mixture containing naphthenic acids (NAs) as the primary acid component [1, 2]. High-resolution Orbitrap mass spectrometer is well suited for the characterization of complex samples and is useful for monitoring the changes in distribution of OSPW-NAs during treatment. In this work, the effectiveness of constructed wetland treatment (CWT) in the transformation of naphthenic acid fraction components (NAFCs) was evaluated in a non-aerated wetland design using high-resolution Orbitrap-MS.

Methods

NAs transformation in non-aerated CWT was assessed using high-resolution Orbitrap-MS technique. The NAs from aqueous wetland outflow samples were extracted using WAX SPE method described previously [3]. Class distributions were established with acquired accurate mass data and Composer version 1.5.3 with mass accuracies of less than 2 ppm. In addition, kinetic and statistical analysis was performed with SigmaPlot software.

Results

The distributions of NAs in the untreated and CWT were examined for selectivity based on the changes in carbon number numbers and double bond equivalents (DBE). The results demonstrated that the non-aerated CWT preferentially transformed O₂-NAs with larger carbon numbers than small carbon number NAs. The oxidized NAFCs (O₃-NAFCs and O₄-NAFCs) were more recalcitrant towards transformation relative to O₂-NAs in the wetland system. Furthermore, kinetic data indicated that rate of transformation of O₂-NAs with larger carbon numbers was faster than small carbon number NA species and there was no direct correlation between the increased transformation rate and size of the DBE observed in the treatment system.

Conclusions

The non-aerated wetland treatment system was effective for treatment of OSPW-NAFCs and favors the transformation of higher carbon number over small carbon number species. This suggests that preferential removal of NA species may be related to design of the system. Work is ongoing to compare with other wetland designs such as aerated wetlands to determine whether these designs have impacts on their selective transformation.

Novel Aspect

The present study reported herein is the first study to evaluate the feasibility of OSPW-NAFCs transformation in constructed wetland treatment (CWT) using high-resolution negative-ion electrospray Orbitrap-MS for detailed understanding of their fates and behavior in the CWT.

References

1. Grewer, D. M., Young, R. F., Whittal, R. M., Fedorak, P. M., *Sci Total Environ*, 408: 5997-6010 (2010).
2. Headley, J. V., Peru, K. M., Mohamed, M. H., Frank, R. A., Martin, J. M., *J Environ Sci Heal A*, 48, 1145-1163 (2013b).
3. Ajaero, C., McMartin, D. W., Peru, K. M., Bailey, J., Haakensen, M., Friesen, V., Martz, R., Hughes, S. A., Brown, C., Chen, H., McKenna, A. M., Corilo, Y. E., Headley, J. V., *Energ Fuel*, 31, 1731-1740 (2017).

ThP-43 / PHOTOCATALYTIC ACTIVITY OF TITANIUM DIOXIDE FILMS UNDER VISIBLE LIGHT EVALUATED BY ESI(-)-MS

Diego Holanda Pereira de Souza - Thiago Carita Correra - Vitor de Moraes Zamarion

Universidade de São Paulo, Institute of Chemistry, São Paulo

Keywords: Photocatalysis; Titania; Visible light; Organic Pollutants; Electrospray.

Introduction

Some metal oxides, such as titanium dioxide (TiO₂), can be employed as photocatalysts because of their ability to promote oxidation of organic molecules when activated by UV-light [1] due to its wide band gap. In this work, we showed via ESI(-)-MS that thin mesoporous films of TiO₂ nanoparticles can degrade hazardous environmental pollutants, as 2-naphtol [2], under visible light incidence.

Methods

Three identical aqueous solutions of 2-naphtol (photosensitive substance) and ammonium hexafluorophosphate (internal standard) were prepared. To one of them, a glass plate was immersed (control), and into the others two plates of TiO₂ films were dived. One of the tubes containing titania was kept in the dark, and the other two tubes were maintained under incidence of blue LED light.

Results

To conduct a quantitative experiment, NH_4PF_6 was employed as an internal standard for the analysis. This compound was chosen because it easily ionizes in aqueous solution generating anions that are stable and inert under these experimental conditions.

The photodegradation rate was carried by periodically sampling a small aliquot of the solutions in contact with the glass plates every 60 minutes for 6 hours. The decay rate of 2-naphtol was controlled by the ratio of the intensities of the $[\text{2-naphtol-H}]^-$ to the PF_6^- over time.

It was observed that TiO_2 films can easily promote photooxidation of 2-naphtol under incidence of visible light and in the presence of oxygen gas by showing an almost complete photodegradation after about 6 hours of reaction. The systems that were set up with TiO_2 film in the dark and with the bare glass under blue light did not demonstrate any kind of degradation.

Conclusions

TiO_2 films can photodegrade organic pollutants, as 2-naphtol, in the absence of UV radiation and under incidence of visible light. A quantitative MS technique was developed and NH_4PF_6 was used as internal standard. This procedure worked well for evaluating photodegradation and allows the degradation products to be detected.

Novel Aspect

The use of ammonium hexafluorophosphate as an internal standard and the determination that visible light can induce photocatalysis in TiO_2 films.

References

1. Ren H. T., Yang Q., Applied Surface Science, 396, 530-538 (2017).
2. Shahrezaei A. R. F., Akhbari A., International Journal Energy and Environment, 3, 267-274 (2012).

ThP-44 / STUDY ON THE DEGRADATION MECHANISM OF ATRAZINE IN SEWAGE BY UHPLC-MS/MS

Jun Li⁽¹⁾ - *Jianfei Yan*⁽¹⁾ - *Bo Lai*⁽¹⁾ - *Xiang Lia*⁽²⁾ - *Kaveh Kahen*⁽³⁾ - *Jingcun Wu*⁽³⁾

Department of Environmental Science and Engineering, School of Architecture and Environment, Sichuan University, Chengdu, China (1) - PerkinElmer Management (Shanghai) Co., Ltd., Shanghai, China (2) - PerkinElmer Inc., Woodbridge, ON (3)

Keywords: degradation mechanism, atrazine removal, nanoparticles, electrochemical reaction, UHPLC-MS/MS

Introduction

CuFe_2O_4 (CFO) magnetic nanoparticles, worked as both particle electrode and catalyst for persulfate (PS), are applied to the three-dimensional electrochemical processes for atrazine removal. In this study, the effects of several key parameters, such as the content of Na_2SO_4 electrolyte, CuFe_2O_4 , persulfate, other inorganic anions, the initial pH and current density on atrazine removal were investigated. In addition, the degradation intermediates were determined by UHPLC-MS/MS. Based on the identified intermediates, the degradation pathway or mechanisms for atrazine in the CFO/PS system were confirmed. Finally, the possible electrochemical reaction mechanism of the CFO/PS system was proposed according to the comprehensive study.

Methods

The degradation intermediates of atrazine were determined by coupling PerkinElmer UHPLC to a QSight™ 220 mass spectrometer. A C18 column was used for separation with gradient elution at a flow rate of 0.3 mL/min, with mobile phases consisted of acetonitrile (A) and pure water containing 0.1% formic acid (B).

The mass spectra with electrospray ionization (ESI) source were recorded across the range of 50-300 m/z with the positive scan mode and MRM mode. The operation parameters were as follows: drying gas of 100, HSID temperature of 250°C, nebulizer gas of 200, electrospray voltage of +5500V, and source temperature of 450°C.

Results

In this electrochemical reaction process, CFO magnetic nanoparticles have been successfully applied as a catalyst for the activation of persulfate, exhibiting excellent performance for atrazine removal in CFO/PS system after 35 min treatment. The optimal operating conditions (*i.e.*, CFO dosage of 3.0 g/L, PS dosage of 4.0 mM, initial pH of 6.3 and treatment time of 35 min) were obtained by investigating the effects of key parameters on the atrazine removal. The high-efficiency of the CFO/PS for atrazine removal was confirmed through seven control experiments. It was found that CFO particles retained high activity after four excessive cycles, demonstrating the stability and recyclability of CFO in the studied system. The comparative study of the freshly prepared and reacted CFO also confirmed the stability of CFO. In this study, the sulfate radical was identified as the main reactive radical in the atrazine degradation process. In addition, the degradation mechanism of atrazine was ascertained according to the 12 degradation intermediates detected by UHPLC-MSMS. Finally, the possible electrochemical reaction mechanism of CFO/PS system was proposed according to the aforementioned analysis.

Conclusions

In this study, CFO magnetic nanoparticles have been successfully applied as a catalyst for the activation of persulfate, exhibiting excellent performance for atrazine removal in CFO/PS system. The high-efficiency of the CFO/PS for atrazine removal was confirmed. It was found that CFO particles retained high activity after four excessive cycles, demonstrating the stability and recyclability of CFO in the studied system. The sulfate radical was identified as the main reactive radical in the atrazine degradation process and the degradation mechanism of atrazine was ascertained according to the 12 degradation intermediates detected by UHPLC-MSMS.

ThP-45 / COMPARISON OF ORGANIC MIXTURES FROM PARTICULATE MATTERS COLLECTED IN KOREA AND CHINA BY USING GCXGC/HIGH RESOLUTION MASS SPECTROMETRY

Moon-Hee Park⁽¹⁾ - **Geul Bang**⁽¹⁾ - **Kyung-Soon Jang**⁽¹⁾ - **Kihong Park**⁽²⁾ - **Sunghwan Kim**⁽³⁾ - **Young Hwan Kim**⁽¹⁾

Korea Basic Science Institute, Biomedical Omics Group, Cheongju (1) - Gwangju Institute of Science and Technology, School of Earth Science and Environmental Engineering, Gwangju (2) - Kyungpook National University, Department of Chemistry, Daegu (3)

Keywords: particulate matter, PM2.5, GCxGC/HRMS, PAHs, atmosphere air

Introduction

In Korea and neighboring China, airborne particulate matter (PM2.5) is very serious environment problems, having primary organic pollutants directly released from emission sources and secondary organic pollutants generated from atmosphere chemical reaction. Thus, the complex organic compounds extracted from PM2.5 collected day after day during one month in each country were analyzed to compare their identifications, relative quantities and emission sources.

Methods

Samples of PM2.5 were simultaneously collected day after day for 28 days (4-31 January 2018) in Gwangju and Beijing using a high volume air sampler. The 90 cm² filter was extracted with dichloromethane through sonication twice. The two days of extracts were combined, filtered and concentrated under N₂ gas.

Comprehensive two-dimensional gas chromatography/high resolution time-of-flight mass spectrometry (LECOGC-HRT) was utilized to analyze the organic extracts. [1]

Results

Approximately, 460 compounds were separated on the polar and sequential nonpolar GC columns and identified based on the mass spectral data from NIST and Wiley libraries, and exact mass accuracy (<1 ppm) of molecular ion from high resolution data, including alkanes, carboxylic acids, hopanes, PAHs, substituted aromatics and steranes, so on. A variety of PAHs were identified in Beijing samples such as 11H-benzo[b]fluorine, benzo[ghi]fluoranthene, benzo[a]anthracene, acepyrene, pyrene and benzo[A]yrene. Their concentrations in Beijing were more 49 to 3 times according to PAH compounds than those in Gwangju for 16/17 days PM2.5. On the contrary to this, natural resin such as dehydroabietic acid was much more in Korea sample. The organic compound markers which can make a distinction between the emission sources of PM2.5 in Korea and China are present. Among them, naphthalic anhydride found in this analysis may become a candidate marker. The study about what this compound in PM2.5 is originated from is currently ongoing. [2]

Conclusions

The complex organic mixtures from PM2.5 were analyzed by using GCxGC/HRMS. The identified organic compounds were different qualitatively and quantitatively between Korea and China. The PM2.5 generated in China has been sometimes transferred into Korea by atmosphere air. Thus, there is anticipated that the issues between two countries related about the emission sources will be discussed, based on these objective results obtained for PM2.5 collected in Korea and China.

Novel Aspect

The organic compounds extracted from PM2.5 collected in Korea and China were different qualitatively and quantitatively due to their different emission sources.

References

1. Alam, M. S., West C. E., Scarlett A. G., Rowland S. J., Harrison R. M., *Atoms. Environ.*, 65, 101(2013).
2. Lee J. Y., Lane D. A., Heo J. B. Yi S. -M., Kim Y. P., *Atom. Environ.*, 55, 17(2012).

ThP-46 / BIOCIDES SPRAY PRODUCT EXPOSURE: MEASURED GAS, PARTICLE AND SURFACE CONCENTRATIONS COMPARED WITH SPRAY MODEL PREDICTIONS

Vivi Kofoed-Sørensen⁽¹⁾ - *Asger W. Nørgaard*⁽²⁾ - *Marie Frederiksen*⁽¹⁾ - *Peder Wolkoff*⁽¹⁾ - *Per Axel Clausen*⁽¹⁾

National Research Center for the Working Environment, Chemical working environment, Copenhagen (1) - Novo Nordisk, Development, Copenhagen (2)

Keywords: Spray products, Biocides, Modelling, Exposure, MS

Introduction

Spray product emissions may be inhaled and some of the chemical ingredients may be respiratory irritants which may cause development of asthma or have other health effects. The aim of this project was to compare experimental exposure of 3 biocidal spray products (#1, #2, #3) applied in an exposure chamber with computed exposure calculations in existing models. The results were used to validate the models to contribute to a better risk assessment of spray products.

Methods

Three biocidal spray products with different active compounds and spray nozzles were selected. They were

applied in a 20.3 m³ controlled exposure chamber in accordance with the description provided by the manufactures. Airborne organic compounds were sampled on Tenax and XAD and surface deposited by wiping with swabs. Samples were analyzed by (TD)-GC-MS/MS and LC-(QTOF)MS. Mass and number concentrations and size distributions of particles were measured on-line.

Results

#1 had permethrin, piperonylbutoxide, and pyrethrum extract as active compounds, which along with alkanes were found in the air after spraying. Like the particles these followed a 1st order decay, however, the removal rates were higher than the air exchange rate (AER) indicating they were associated with the depositing particles. The alkane decay followed the AER indicating they were in the gas-phase.

#2 contained λ-cyhalothrin as the active compounds but it was not found in the air samples. Large particles (>0.5 μm) were found in very low concentrations (<0.1 mg/m³), while smaller particles (<0.3 μm) stayed at background levels. λ-cyhalothrin was found in wipe samples from surfaces.

#3 had benzalkonium chlorides as the active compounds, which were measureable in air up to 36 min after spraying following a 1st order decay. The particle pattern was the same as for #2. The wipe samples showed that the benzalkonium chlorides were mainly on surfaces near the spray area.

Preliminary results showed discrepancy between models and experiments.

Conclusions

The sprayed aerosols had shorter residence time in air than explained by the AER. Concentrations of the active compounds and low vapor pressure organic compounds decayed faster than explained by AER since they were associated with the aerosols, while more volatile compounds followed the AER. Aerosols of benzalkonium chlorides could be measured in the air more than ½ hour after application. The active compounds deposited mainly at the floor near the spraying position.

Novel Aspect

Modelling is used extensively in risk assessment; however there is a research gap on the compliance between modelling and actual exposure data using standardized conditions.

ThP-47 / OCCURRENCE AND HUMAN EXPOSURE TO BROMINATED AND ORGANOPHOSPHATE FLAME RETARDANTS FROM FINNISH HOUSEHOLD DUST

Eva Kumar⁽¹⁾ - *Panu Rantakokko*⁽¹⁾ - *Joris Braber*⁽¹⁾ - *Jani Koponen*⁽¹⁾ - *Taya Huang*⁽¹⁾ - *Hannu Kiviranta*⁽¹⁾

National Institute for Health and Welfare (THL), Health Security, Kuopio, Finland (1)

Keywords: brominated flame retardants, organophosphate flame retardants, household dust, human exposure, exposure model

Introduction

Flame retardants (FRs) are chemicals that are used to decrease flammability of consumer products. They are incorporated as additives or reactive ingredients in materials including plastics and textiles. Brominated (BFRs) and organophosphate flame retardants (PFRs) are the two main classes of FRs that are extensively used wherein the latter are considered as potential replacements to conventional BFRs but also part of those exhibit potential health risk to consumers.

Methods

BFRs and PFRs were measured in Finnish household dust (n=40) by GC-HRMS. Children and adults exposure assessment (ingestion, inhalation, dermal) was based on measured dust concentrations, and it was made according to models published by Little et al. [1]. Exposure assessment needed calculation of dust-air distribution of FRs, the validity of which was evaluated using measured dust and indoor air concentrations of FRs from Norway [2].

Results

The certified dust SRM2585 was used as the quality control of the analytical method: recoveries were 91-121% for BFRs and 65-120% for PFRs. In case of BFRs, BDE-209 (410 ng/g) had the highest median concentration followed by DBDPE (Decabromodiphenyl ethane) (120 ng/g), and BEH-TEBP (Bis(2-ethylhexyl)tetrabromophthalate) (110 ng/g) in household dust. TBOEP (Tris(2-butoxyethyl)phosphate) (11000 ng/g) was the most abundant PFR followed by TCIPP (Tris(2-chloroisopropyl)phosphate) (1870 ng/g) and TPHP (Triphenyl phosphate) (770 ng/g). From modelled exposure assessment, it was determined that exposure via dust ingestion contributes most for BFRs and also for less-volatile PFRs for children who are also more susceptible to FRs exposure than adults. But dermal adsorption governs the exposure pathway for more-volatile PFRs. However, the model used for determining air concentrations under-estimated some low volatile BFRs and over-estimated some volatile PFRs, that results in the respective under/over estimation of exposure.

Conclusions

Measured dust and calculated air concentrations of FRs in Finland were found to be similar with those measured in Norway. Total daily FR intakes were some orders of magnitude lower than the oral reference doses from literature. This study presents simple analytical and exposure assessment methods which enable effective screening of exposure. However, further studies are required to develop cost-effective methods to study the distribution and exposure of FRs.

Novel Aspect

We present a simple analytical-calculation method to assess total exposure from indoor environment to FRs that is based solely on the measured FR concentrations in dust.

References

1. Little J. C., Weschler C. J., Nazaroff W. W., Liu Z., Hubal E. A. C., Environ Sci Tech, 46, 11171 (2012).
2. Cequier E., Ionas A. C., Covaci A., Marcé R. M., Becher G., Thomsen C., Environ Sci Tech, 48, 6827 (2014).

ThP-48 / NOVEL DRINKING WATER DISINFECTION BY-PRODUCTS

Albert Lebedev^(1,2) - ***Dmitry Kosyakov***⁽²⁾ - ***Nikolay Ul'yanovskii***⁽²⁾ - ***Mark Popov***⁽²⁾ - ***Tomas Latkin***⁽²⁾

Moscow State University, Moscow, Russian Federation (1) - Northern (Arctic) Federal University, Arkhangelsk, Russia (2)

Keywords: disinfection by-products, HPLC-HRMS, tap water, organochlorines, organobromines

Introduction

Drinking water chlorination is used over 100 years. Destroying microbes chlorination however leads to formation of toxic disinfection by-products (DBPs), appearing due to reaction of active chlorine with natural and anthropogenic compounds dissolved in water. The list of DBP, contains over 700 compounds. However they cover only 50% of organic halogens in drinking water. In the present study over 30 new halogenated organic compounds were identified by HPLC-HRMS.

Methods

Water samples from the Severnaya Dvina river at the water intake of Arkhangelsk (Russia), after purification and coagulation, and after chlorination were collected in March-June 2017. Two liters of water at pH 3 were made pass through SPE Bond Elute PPL cartridges (Agilent) cartridges. High resolution QTOF mass spectrometer TripleTOF 5600, (AB Sciex) with ion source Duospray connected to the HPLC system LC-30 Nexera (Shimadzu) was used.

Results

Non targeted screening of pollutants in the Arkhangelsk tap water with HPLC-HRMS allowed detecting over 50 organohalogenes. Among phenols pentachlorophenol is usually considered as an independent priority pollutant. However it was not present in the samples before chlorination.

MS/MS allowed referring the largest group to halohydrines and dihalogenated fatty amides, the products of conjugated electrophilic addition of halogens to unsaturated fatty amides. The proposed structures were confirmed by oleamide chlorination in water. These compounds are a brand new class of DBPs, while their toxicities require special study. Five oxygenated species were 3,5-dichloro-4-hydroxybenzoic, dimethyl-bromobenzoic, and dimethyldibromobenzoic acids, trichlorophenylacetate, and 2,4-trichloro-1,3,5-benzotrialdehyde. ChemSpider search and spectra interpretation of *N*-containing compounds allowed identifying 3,4-dichloro-2-methoxy-5-oxo-2,5-dihydro-1H-pyrrol-1-yl)-acetic acid and highly substituted halogenated pyridines.

Conclusions

Several groups of new disinfection by-products were identified in the tap water of Arkhangelsk (Russia). Mass spectral conclusions were confirmed by independent synthesis with standard compounds. Chlorinated and brominated compounds are formed at the stage of disinfection of the Severnaya Dvina river water rich with humic matter. Unfortunately there is no any available information on the toxicities of the novel DBPs.

Novel Aspect

Several groups of novel disinfection by-products were established with HPLC-HRMS in the drinking water in Arkhangelsk, Russia

ThP-49 / A COMPARATIVE STUDY OF ULTRA-TRACE-LEVEL URANIUM BY THERMAL IONIZATION MASS SPECTROMETRY WITH CONTINUOUS HEATING: STATIC AND PEAK-JUMPING MODES

Chi Gyu Lee - Jinkyu Park - Sang Ho Lim - Sun Ho Han

Korea Atomic Energy Research Institute, Nuclear Chemistry Research Division, Daejeon

Keywords: Ultra-trace-level Uranium, Isotope Ratio, TIMS, Detection mode, Nuclear Safeguards

Introduction

For ensuring nuclear safeguards, we report the analytical signal-detection performance of thermal ionization mass spectrometry (TIMS) with continuous heating for the measurement of isotopic ratios in samples containing ultra-trace amounts of uranium.

Methods

The isotope ratio measurements were performed using a thermal ionization mass spectrometer (TRITON Plus, Thermo Fisher Scientific, Germany). As methods for detecting uranium signals, peak-jumping mode using a single detector and static mode using multiple detectors were examined with U100 (10% ²³⁵U-enriched) uranium standard samples in the femtogram-to-picogram range.

Results

Uranium isotope ratios, $n(^{235}\text{U})/n(^{238}\text{U})$, were measured down to levels of 1 fg and 3 fg in static and peak-jumping modes, respectively, while $n(^{234}\text{U})/n(^{238}\text{U})$ and $n(^{236}\text{U})/n(^{238}\text{U})$ values were measured down levels of 100 fg in both modes. In addition, the dependency of the ²³⁸U signal intensity on sample quantity exhibited similar tendencies in both modes. The precisions of the isotope ratios obtained in the static mode over all sample ranges used in this study were overall slightly higher than those obtained in peak-jumping mode.

Conclusions

These results indicate that the analytical performance by TIMS with continuous heating are almost independent of the detection method, *i.e.*, peak-jumping mode or static mode, which is characteristic of isotope-ratio measurements using the TIMS method with continuous heating.

Novel Aspect

The results suggest that similar analytical performance can be obtained in the two detection of TIMS with continuous heating.

ThP-50 / ANALYSIS OF NATURAL RADIONUCLIDES CONCENTRATIONS IN MONAZITE USING ACID DIGESTION BY QUADRUPOLE ICP-MS

Juhyun H. Lee ^(1,2) - *Sung Yi* ⁽²⁾

Korea Institute of Nuclear Safety, Department of Natural Radiation Safety, Daejeon, Republic of Korea (1) - Hanyang University, Department of Chemical Engineering, Seoul, Republic of Korea (2)

Keywords: Acid digestion, uranium, thorium, ICP-QMS

Introduction

According to the Act on Protective Action Guideline Against Radiation in the Natural Environment in South Korea has been planned and conducted the field investigation to monitor the radioactivity concentration surrounding usage facilities of raw materials and residues and protect from radiation in the natural environment. Anyone who handles raw materials or residues contained uranium or thorium whose activity concentration exceeds 1 Bq/g should be registered the material informations (radioactivity concentrations, total amount, etc.). The objective was to monitor of potential pollution materials in naturally environment and offer some useful information for policy decision. Then it was important to pretreat the samples in order to accurately analyze the concentrations. There are various digestion procedures: open vessel digestions, closed vessel digestions, microwave digestions, alkali fusion. The digestion of monazite is an important aspect of thorium and uranium elements analysis using by quadrupole inductively coupled plasma mass spectrometry (ICP-MS). This paper describes comparison of two digestion methods for thorium and uranium elements analysis in monazite.

Methods

In this study, the monazite was dissolved by the conventional acid digestion and electric fusion system. Prior to the pretreatment, it was using a XRD in order to analyze the components of the suspicious materials, detected matter whose radioactive concentration exceeds or is suspected to exceed limits determined by NSSC in monitors at harbors. After homogenous sample preparation through milling or mixing, the sample were dried in a temperature controlled oven up to 450°C for about 5hrs. The samples were dissolved by fusion with lithium metaborate and acid digestion with HF, HNO₃, and HClO₄ to break refractory complex. The dissolved solutions were precipitated at pH 7 using Fe co-precipitation and dilution HNO₃ for analysis by ICP-QMS. In addition, the samples were completely sealed in D6H4 bottle to secure the confinement of radon gases. It spends time for 1 month to attain the secular equilibrium between radium and radon isotopes in decay series. Thorium and uranium series that represented 228Ac and 214Bi gamma nuclides measured by HPGe and ²³²Th and ²³⁸U atomic analyzed by ICP-QMS.

Results

The reference sample with certificate values (zircon by ALML) was analyzed for the different digestions. The alkali fusion method showed great thorium and uranium yield up to 98% and the closed vessel acid digestion method was similar thorium and uranium yield. The results of ²³²Th and ²³⁸U using HPGe and ICP-QMS showed a similar distribution. However, the monazite showed different values between fusion

method and acid digestion. The higher concentration of the samples, the bigger the difference values. In addition, the radioactivity of the ^{232}Th and ^{238}U indicator decreased in the monazite and showed a different concentration of thorium and uranium between ICP-QMS and HPGe.

Conclusions

It is the most important to monitor environmental radiation, raw materials, residues and products for protect from radiation exposure. NSSC enforced "Act on Protective Action Guidelines Against Radiation in the Natural Environment" since 2012. Thus, this study establishes the procedures of radiation measurement or radioactivity analysis about NORM samples. In particular, it was necessary to analyze the thorium and uranium concentrations of monazite using acid digestion and evaluate the assessment of radiation risk for the Korean population and protect the health of the workers involved and the environment.

Novel Aspect

Comparison between alkali fusion and acid digestion for thorium and uranium analysis in monazite or products suspected as contaminant materials in the field survey based on NORM regulation in South Korea

References

1. Y. G. Ha., Journal of the Korean Chemical Society, vol. 23, No. 3 (1979).
2. L. Kotz, G. Kaiser, P. Tschöpel und G. Tölg Z. Anal. Chem. 260, 207-209 (1972).

ThP-51 / DETERMINATION OF GLOBAL DNA METHYLATION AND OXIDATION OF DNA AND LIPID BIOMARKERS IN METAL OXIDE NANOMATERIAL-HANDLING WORKERS BY LIQUID CHROMATOGRAPHY-TANDEM MASS SPECTROMETRY

Hui-Ling Lee⁽¹⁾ - *Saou-Hsing Liou*⁽²⁾

Fu Jen Catholic University, Fu Jen Catholic University/Chemistry, New Taipei City (1) - National Health Research Institutes, National Institute of Environmental Health Sciences, Zhunan, Miaoli County (2)

Keywords: metal oxide nanomaterials; 8-hydroxydeoxyguanosine; biomarkers; LC-MS/MS

Introduction

The nanotechnology industry is a rapidly growing. The most commonly used nanomaterials (NMs) are, titanium dioxide (TiO_2), silica (SiO_2) and zinc oxide. To assess whether oxidative stress markers of nucleic acids (8-OHdG), DNA methylation levels in urine and white blood cell (WBC), and EBC 8-isoprostane were associated with nanomaterial exposure. Moreover, the relationship between oxidative damage biomarkers and global DNA methylation was investigated.

Methods

In the present study, urinary and white blood cell (WBC) 8-hydroxydeoxyguanosine (8-OHdG), and exhaled breath condensate (EBC) 8-isoprostane were measured as oxidative stress biomarkers. WBC global methylation was measured as an epigenetic alteration. To examine global methylation, oxidative DNA damage, and lipid peroxidation in workers with occupational exposure to metal oxide nanomaterials using LC-MS/MS.

Results

Exposure to TiO_2 , SiO_2 , and indium tin oxide (ITO) resulted in significantly higher oxidative biomarkers such as urinary 8-OHdG and EBC 8-isoprostane. However, significantly higher WBC 8-OHdG and lower global methylation were only observed in ITO handling workers. Significant positive correlations were noted between WBC and urinary 8-OHdG (Spearman correlation $r = 0.256$, $p = 0.003$). Furthermore, a significant

negative correlation was found between WBC 8-OHdG and global methylation ($r = -0.272$, $p = 0.002$). These results suggest that exposure to metal oxide NMs may lead to global methylation, DNA oxidative damage, and lipid peroxidation.

Conclusions

Global DNA methylation levels, and levels of oxidative stress markers 8-OHdG and 8-isoprostane were assessed in metal oxide nanomaterial handling workers. The results show that exposure to metal oxide nanomaterials may lead to global hypomethylation, DNA oxidative damages, and lipid peroxidation. Furthermore, we identified a negative correlation between WBC 8-OHdG and global methylation.

Novel Aspect

Exposure to metal oxide nanoparticles may lead to global methylation and DNA oxidative damage and lipid peroxidation.

References

1. Wang C. J., Yang N. H., Chang C. C., Liou S. H., Lee H. L., *Journal of Chromatography B* 879, 3538-3543 (2011).
2. Wang C. J., Yang N. H., Liou S. H., Lee H. L., *Talanta* 82, 1434-1438 (2010).

ThP-52 / MASS SPECTROMETRY IDENTIFICATION OF HALOGENATED PEPTIDES IN TREATED DRINKING WATER

Ping Jiang - Guang Huang - Lindsay Blackstock - Xing-Fang Li

University of Alberta, Blackstock, Edmonton, Canada

Introduction

Ionic liquids (ILs) are increasingly used in industrial processes as “green chemicals” because of unique properties of low volatility and customizability. ILs can be used to enable novel processes and/or replace conventional organic solvents in a wide variety of applications. Widespread use may increase the risk of accidental release of IL-containing industrial wastes into environmental waters. Most ILs are highly water soluble, and have estimated environmental half-lives of several days to a month. IL cations often consist of aromatic or alkyl quaternary amines that resemble previously confirmed *N*-nitrosamine (NAs) precursors. NAs are confirmed animal carcinogens and classified as probable human carcinogens. NAs are also potent, estimated to have negative health effects at ng/L concentrations. Therefore, this study is to investigate the formation of ILs of various structures to form nitrosamines under the conditions of water chloramination.

Methods

Laboratory reactions of ILs consisting alkyl and aromatic tertiary amines with chloramine will be performed. HPLC-MS and HPLC-MS/MS methods are developed and used to identify and quantify the nitrosamines produced in the reactions.

Preliminary Results

We have studied two ILs, 1-ethyl-1-methylpyrrolidinium bromide (EMPy) and 1-ethyl-3-methylimidazolium bromide (EMIm), for their nitrosamine formation potential. Each IL species was reacted with pre-formed monochloramine under various conditions. After 24h, samples were extracted from water using dichloromethane. The extracts were analyzed for nitrosamines formed in the reactions. HPLC-high resolution MS analysis of the extracts identified the produced nitrosamines, while HPLC-MS/MS with multiple reaction monitoring mode quantified the concentrations of the NAs using a deuterated internal standard. Both EMIm and EMPy can produce NAs: EMIm generating *N*-nitrosomethylethylamine, while EMPy producing *N*-nitrosopyrrolidine. The formation yield of NAs is dependent on the structures of the ILs.

EMPy_r produced higher amount of NA than EMIm under all conditions evaluated. EMPy_r can produce NA with a yield on the same order of magnitude as polydiallyldimethylammonium chloride, a confirmed nitrosodimethylamine precursor. We are investigating a variety of ILs for the formation potential of NAs. In addition, we will further investigate how to use UV/H₂O₂ to reduce the formation of nitrosamines from ILs. This study is the first report demonstrating the formation of NAs from ILs under water chloramination conditions. This study emphasizes the importance of prevention of environmental discharge of ILs to water bodies, and highlights a need for further evaluation of potential lifecycle impacts of ILs prior to their wide ranging applications.

Novel Aspect

MS studies of Ionic liquids forming toxic nitrosamines under common water treatment conditions.

ThP-53 / IODINE DISTRIBUTION AND BIOAVAILABILITY IN SOIL FROM FENYANG, SHANXI PROVINCE, CHINA

Wei Liu - Jundong Hu

Chinese Academy of Geological Sciences, National Research Center for Geoanalysis, Beijing

Keywords: iodine, speciation, stability, distribution, soil

Introduction

Soil plays an important role in the environmental geochemistry circulation and transformation process of iodine. In this study, the primary interest was to understand the iodine species that were potentially mobilizable and bioavailable, rather than the total iodine species present in the soil. Investigation has been conducted on the status of iodine distribution in Fenyang, in order to elucidate the distribution characteristics of iodine in the environment.

Methods

Soil profile samples were mainly collected from farmland in Fenyang. The total iodine concentration was detected by inductively coupled plasma mass spectrometry (ICP-MS) using ammonium hydroxide digestion. The speciation analysis of iodine was determined by high performance liquid chromatography (HPLC) coupled with ICP-MS. Two methods that water extracting and phosphate extracting were used to investigate the bioavailability of iodine in soil samples.

Results

Good results were achieved using 10 % ammonium hydroxide (v/v) and 0.5 mM EDTA where the transformation between IO₃⁻ and I⁻ was clearly minimized. Inorganic iodine was the main iodine species in the extraction solution. Iodine in soil exists as organic iodine has been reported elsewhere [1,2,3]. Our results indicated that the transformation of organically bound iodine to inorganic iodine must be occurred because the sums of inorganic iodine were agreed with the certified values of reference soils. But the transformation between IO₃⁻ and I⁻ was not observed.

Water extraction and phosphate extraction were used to investigate the bioavailability of iodine in 6 soil samples. The results showed that I⁻ was the main iodine species. It confirmed that iodine in soils was mainly transported into plants in the form of inorganic iodine. The iodine concentrations of all the soil profile samples ranged from 0.56~5.69 µg-I g⁻¹. The maximum iodine concentration was at the depth of 0~20 cm. And the soil iodine content decreased with increasing sampling depth.

Conclusions

Efforts were made to preserve the oxidation of IO₃⁻ to I⁻ by using EDTA. Although this method should be further tested for its application for a broad range of soil samples, the results obtained in this study are promising. The main species of iodine uptake from soil was iodide. The content of iodine in the surface soil

was higher than that in the depth. It was not only related to iodine content in groundwater, but also related to geological structure and properties of soil.

Novel Aspect

The stability of inorganic iodine was achieved by adding EDTA into the extraction solution. The distribution and bioavailability of iodine in soil profile samples was discussed.

References

1. Schlegel, M. L., Reiller, P., Bion, F. M., Barré, N., Moulin V. *Geochimica et Cosmochimica Acta*. 70, 5536-5551 (2006).
2. Shimamoto, Y. S., Takahashi, Y., Terada, Y. *Environmental Science and Technology*. 45, 2086-2092 (2011).
3. Yamaguchi. N., Nakano, M., Takamatsu, R., Tanida, H. *Journal of Environmental Radioactivity*. 101, 451-457 (2010).

ThP-54 / DETERMINATION OF Ag NANOPARTICLES IN VLTAVA RIVER (PRAGUE, CZECH REPUBLIC) BY SINGLE PARTICLE ICP-MS

Martin Loula - Antonin Kana - Oto Mestek

University of Chemistry and Technology Prague, Department of Analytical Chemistry, Prague

Keywords: single particle ICP-MS; silver nanoparticles; Vltava river

Introduction

Nanoparticles are used in many areas of industry and everyday life and they are released in great quantities into the environment. Despite this, little is known about their (eco)toxicity and the risks they might pose. To be able to evaluate these risks and take appropriate precautions, it is necessary to understand their behavior and fate in the environment.

Methods

The single particle ICP-MS was used to determine and characterize Ag nanoparticles in samples of river water taken at several places from Vltava and a few smaller streams in Prague and its neighborhood. NexION 350D (Perkin Elmer, Concord, Canada) spectrometer was used for all measurements. Samples were stabilized by addition of gelatin immediately after sampling and were analyzed the same day.

Results

The developed method exhibited good characteristics and only minor matrix effects of the river water were observed: limit of detection (LOD) of particle size ca. 25 nm (depends on dissolved Ag concentration), LOD of particle number concentration 100 mL⁻¹, 1.3% repeatability and 98% recovery. Small distortion of nanoparticle diameter was observed: mean diameter 92 nm was found for particles with diameter 96 ± 8 nm.

Ag nanoparticles were present in samples of water from Vltava river. The found Ag nanoparticles covered a wide range of diameters, ca. 30-120 nm. Their number concentrations were extremely low, ca. 300-2000 mL⁻¹, however well above the LOD. The Ag nanoparticles number concentrations increased as the river flowed through the urban agglomeration.

Conclusions

Single particle ICP-MS proved to be capable of determining Ag nanoparticles in surface waters. Ag nanoparticles with wide range of diameters were successfully determined in the Vltava river and an increase of number concentrations, likely due to human activity, was observed downstream of Prague. It should be evaluated whether these number concentrations can be harmful to the environment.

Novel Aspect

To the best of our knowledge this is the first study of Ag nanoparticles in surface waters in the Czech Republic.

Acknowledgements

Financial support from specific university research (MSMT No 21-SVV/2018) and Grant Agency of Czech Republic (project no. 17-00291S) is gratefully acknowledged.

ThP-55 / ANTIBIOTICS IN WATERS ALONG THE GRAND CANAL, CHINA

Guohui Lu - Haitao Piao - Shu Chen - Xingchun Jiao - Yongliang Yang

Chinese Academy of Geological Sciences, National Research center for Geoanalysis, Beijing

Keywords: Antibiotics; river water; the Grand Canal; spatial distribution

Introduction

The Grand Canal passes through some industrial cities (e.g. Hangzhou, Suzhou, Wuxi, and Tianjin), and also vast rural areas. The aim of this work is to reveal the differences in contamination levels and composition profiles of antibiotics in waters of different sections of the Canal with comparisons of population density and economic scale in each section.

Methods

HLB column (500 mg, 6 cc, Waters, USA) was used as the solid phase extraction for extraction and enrichment of the target antibiotics. Analysis of antibiotics was performed using a high-performance liquid chromatograph tandem mass spectrometer. Separation of the analytes was performed by an Agilent HP1200 liquid chromatograph interfaced with API 4000 tandem quadrupole mass spectrometer operated in the electrospray ionization mode.

Results

Among the 20 target antibiotics, 14 antibiotics were detected. The predominant types of antibiotics were sulfonamides and macrolides. The spatial distribution pattern of antibiotics shows that the concentration levels at the sampling sites in the north end of the canal (average: 117.1 ng/L) and the south end sections (average: 96.5 ng/L) of the canal were obviously higher than those in the middle sections (average: 48.8 ng/L). At the north end of the canal there are two metropolitan cities Beijing and Tianjin, and at the south end of the canal there are cities with dense population and industries such as Suzhou and Hangzhou. Municipal sewage could be the main sources for sulfonamides and macrolides in these cities[1-2]. Most of the sampling sites in the middle part of the canal are in rural areas, especially in the northern Jiangsu Province, and Shandong and Hebei Provinces where animal husbandry and fish farming could be the main sources for sulfonamides and macrolides.

Conclusions

High concentrations of sulfadiazine, sulfamethoxazole and ampicillin were observed. Antibiotics showed large geographical variations from the North to the South, with high concentrations occurring at the sampling sites in the Yangtze River Delta area and Tianjin-Beijing urban and suburban areas, while low levels were generally found in the middle sections of the canal with vast agricultural areas, showing a V-shape geographic distribution along the river.

Novel Aspect

The impacts of urbanization and economic development on concentrations and compositions of antibiotics in waters of different sections of the Grand Canal are discussed.

References

1. Gao, J., Huang, J., Chen, W., Wang, B., Wang, Y., Deng, S., Yu, G., *Front. Environ. Sci. Eng.* 10,491–501(2016).
2. Dai, G., Wang, B., Fu, C., Dong, R., Huang, J., Deng, S., Wang, Y., Yu, G., *Environ Sci Process Impacts.* 18, 445–455

ThP-56 / GC×GC-MS IN ENVIRONMENTAL ANALYSIS: FROM EFFICIENT SEPARATION TO INTRODUCTION OF RETENTION INDEX DATA BASE

Dmitrii Mazur^(1,2) - *Igor Zenkevich*⁽³⁾ - *Viatcheslav Artaev*⁽⁴⁾ - *Olga Polyakova*⁽¹⁾ - *Albert Lebedev*⁽¹⁾

Lomonosov Moscow State University, Organic Chemistry Department, Moscow (1) - Lomonosov Northern (Arctic) Federal University, Core Facility Center "Arktika", Moscow, Arkhangelsk (2) - St. Petersburg University, Institute for Chemistry, St. Petersburg (3) - LECO Corporation, St. Joseph, MI (4)

Keywords: GC×GC-MS, environmental mass spectrometry, retention index, semivolatile compounds

Introduction

Gas chromatography-mass spectrometry (GC-MS) is one of the most precise, well developed, and reliable analytical tools for the analysis of small volatile and semivolatile compounds. This method has been extensively improved and tested by enhancing the separation technique via comprehensive two-dimensional gas chromatography (GC×GC-MS). Application of GC×GC-MS is especially useful for analysis of complex mixtures, which environmental objects usually are.

Methods

Various organic compounds as well as dichloromethane extracts of different water samples were analyzed using Pegasus® GC×GC HRT 4D instrument (LECO Corp., USA). Compounds from SV Internal Standard Mix (Restek, USA) were used as internal standards (IS) for retention index calculations (perdeuterated 1,4-dichlorobenzene and some PAH). Chromatographic separation was performed using Rxi-5MS and Rxi-17SilMS columns (Restek, USA) for both dimensions respectively.

Results

Several examples revealing insufficiency of GC-MS and GC-HRMS methods in precise identification of some organic compounds put forward GC×GC-MS as one of the most powerful tools for non-target analysis. To improve the identification certainty a simple method for calculation of retention indices attributed to the secondary dimension column (2I) using single injection and PAH as internal standards commonly used in US EPA methods for the analysis of semivolatile organic compounds was proposed. Using the linear combination of both retention time (tR) and retention temperature (TR): $2I = a(TR+ktR) + b$, a series of 2I values for various organic compounds were calculated with good reproducibility under different separation conditions. a and b are proportional to enthalpy and entropy of interaction between analytes and stationary phases, k – non-linearity compensation coefficient. Time constituent of the 2I function by itself was found useful for identification of homological series of some classes of compounds (alkanes, fatty acids).

Conclusions

GC×GC-MS was demonstrated to be crucial for non-target analysis of complex environmental mixtures. A simple and reliable approach is proposed for the calculation of retention indices in GC×GC experiments. The proposed method uses the Lee indices based on PAH, that may be used as a reference under certain conditions, and doesn't require any additional runs. The method is applicable for both target and non-target analyses for a wide range of compounds.

Novel Aspect

A simple and reliable approach is proposed for the calculation of retention indices in GC×GC analyses as one more analytical parameter to increase the identification certainty.

ThP-57 / THE ANALYSIS OF ENVIRONMENTAL SAMPLES USING HIGH RESOLUTION MASS SPECTROMETRY TO IDENTIFY NOVEL PFAS COMPOUNDS

Daniel McMillan⁽¹⁾ - Simon Roberts⁽²⁾ - Craig Butt⁽²⁾ - Chris Borton⁽²⁾

SCIEX, Business Development, Warrington (1) - SCIEX, Application Development, Framingham (2)

Keywords: PFAS, EPA, Non-target, Identification, Q-TOF

Introduction

There are a very wide number of PFAS compounds that are actually present in environmental samples that are not part of the current EPA methods for PFAS analysis. The goal of this work was to use a QTOF mass spectrometer with Non Target Screening workflow to locate and identify unique PFAS compounds that were not known to be in the sample.

Methods

Water samples were analyzed using a Sciex X500R QTOF system set up to acquire data in a Non Target Screening workflow. The instrument first performed a TOF scan from 100-1500. This scan was followed by a data dependent scan where ions above a set threshold were selected for fragmentation and acquisition of a full scan MS/MS spectrum. The resulting MS/MS spectra were searched against a high resolution library to identify potential PFAS of compounds present in the samples.

Results

Perfluoroalkyl substances (PFASs) encompass a range of fully fluorinated alkyl compounds and are prevalent in Aqueous Film Fire Foam (AFFF). In addition, PFASs are ubiquitous as they are used in many household goods and have been found in various environmental and biological samples. Here, we demonstrate the use of QTOF technology to exploit the power of high resolution mass spectrometry and the use of product ion spectra to identify novel compounds in AFFF contaminated samples.

Novel Aspect

Novel PFAS compounds are identified using high resolution mass spectrometry

ThP-58 / CHARACTERIZATION OF ORGANIC COMPOUNDS IN AQUATIC SEDIMENT BY ION PROFILES OBTAINED USING GC/MS SYSTEM

Zaharie Moldovan

National Institute of Research and Development for Isotopic and Molecular Technology, Mass Spectrometry, Cluj-Napoca, Str. Donat 67-103, Romania

Keywords: GC/MS, river sediments, ion profiles

Introduction

The analysis of sediments is an efficient way of re-constructing the depositional history of organic compounds and environmental changes in aquatic environment [1-3]. The organic compounds are released to the environment through anthropogenic (industrial releases) and natural sources (aquatic organisms and

terrestrial plants). Some compounds in sediments are widely used as a sensitive indicator for source identification. In respect to this, one of priority tasks is to determine the environmental distribution of organic compounds in the aquatic sediments. The present paper purpose is the characterization of main families of compounds by characteristic ion obtained by Ion Impact Mass Spectrometry using GC/MS system.

Methods

Sediment sample were prepared by Ultrasonic Solvent Extraction (USE). An aliquot (2g) of freeze-dried sediment was extracted successively with 10 ml CH₂Cl₂+methanol (2:1). The extract was fractioned on Al₂O₃ with following organic phases: n-hexane, n-hexane+CH₂Cl₂ (1:2), CH₂Cl₂+methanol (1:1) and methanol respectively. The final extracted was dried to 1 ml in n-hexane. The analyses were performed using a GC/MS system with mass spectrometer operated in EI mode to 70 eV. The gas chromatograph was equipped with a capillary column HP-5MS (30x0.25mm) with 0.25 film thickness with programmed temperature from 90°C to 315°C.

Results

The families of compounds are visualized by chromatograms on diagnostic ions. The main group of compounds is: 1) *n*-Alkanes which show compounds with the number of carbon atom ranging from 12 to 30.2) Acyclic isoprenoids compounds with a large number of methyl branching points in their structure. This group of compounds include pristine (C19) and phytane (C20) often used as an indicator of depositional environments. 3) Pentacyclic triterpanes (Hopanes) a class of compounds that has a great impact on petroleum geochemistry, they being useful as source indicators. 4) An important group of compounds found in the sediment samples is the Lineal Alkylbenzens (LABs), with number of carbon from 10 to 13.5) An abundant class of compounds consists from Polyaromatic Hydrocarbons (PAHs) of environmental concern due to their carcinogenic proprieties. 8) Sulphonated Polyaromatic Hydrocarbons (SPAHS), other important group of pollutants used very often as biomarkers for the source information.

Conclusions

Are obtained characteristic profiles of organic compound families present in sediments. Their profile shapes is a function of their sources and reflect the natural, industrial and domestic activities at regional level.

Novel Aspect

An efficient method based on GC/MS technique is developed for detection of eight families of organic compounds in sediment. Every family is presented as characteristic ion profile having a large number of compounds.

References

1. Batagliion G. A. et al., Organic Geochemistry 92, 24-31 (2016)
2. Matic I. et al. J. Chromatogr. A 1364, 117-127 (2014)
3. Wang S. et al., Ecotoxicology and Environmental Safety 150, 199-206 (2018)

ThP-59 / STUDY OF PHARMACEUTICALS REMOVAL RATES IN WASTEWATER TREATMENT PLANT USING GC/MS ANALYTICAL TECHNIQUE

Zaharie Moldovan - Olivian Marincas - Ioana Feher

National Institute of Research and Development for Isotopic and Molecular Technology, INCDTIM, Cluj-Napoca, Str. Donat 67-103, Romania

Keywords: GC/MS, PPCPs, WWTP, isotopic dilution, removal rate,

Introduction

Many xenobiotics including Pharmaceuticals and Personal Care Products (PPCPs) are released into the environment from Wastewater Treatment Plants (WWTPs) [1]. The pharmaceuticals persist in the environment, remaining active and affecting aquatic life [2-4]. The PPCPs could act as persistent compounds, because of their continual discharge into aquatic sustaining a multigenerational exposure for the resident organisms [3]. In recent years, reliable methods have been established for analysis of these pollutants down to low ng/L levels in water environment, based on HPLC-MS and GC-MS. The purpose of this work was to quantify some pharmaceutical compounds in the aqueous phase at different stages of a conventional WWTP situated in Cluj-Napoca (Romania) with aim to determination of removal rate.

Methods

The compounds were extracted by SPE procedures and derivatised with *N*-methyl-*N*-trimethylsilyltrifluoroacetamide (MSTFA) and than analyzed by Gas Chromatography-Electron Ionization Ion Trap Mass Spectrometry (GC/EI-ITMS), in full scan mode. For quantification isotopic dilution was applied.

Results

An effective analytical method for the simultaneous determination of pharmaceuticals from various therapeutic classes in aqueous samples has been developed. The optimized GC separation was obtained for 38 PPCPs. The complete mass spectra were registered for natural, derivatised and isotopic labeled compounds. Are detected also compounds originating from thermal decomposition of pharmaceuticals during GC separation. Were analyzed sample collected from different places of WWTP leading to removal rate for different treatment units.

Conclusions

Were obtained: 1) GC/MS separation of 38 compounds from family PPCPs; 2) Mass spectra for natural, silylated and isotopic labeled compounds; 3) Mechanism of thermal decomposition for cyclophosphamide and carbamazepine; 4) The removal rate for every compound in WWTP in the range 19.5-92 %, depending of chemical structure.

Novel Aspect

a) Detection and quantification of PPCPs in an efficient analysis by SPE extraction, derivation by MSTFA and GC/MS analyses; b) Registration of mass spectra of the natural, silylated and isotopic labeled for the same compounds; c) Determination of the removal rate for different units of WWTP.

References

1. Buchberger W. W, J Chromatogr. A, 1218 (2011) 603–618;
2. Fedorova G. et al., Chemosphere 111 (2014) 55–60;
3. Daughton C. G and Ternes T. A. (1999), Environ Health Persp 907-938;
4. Afonso-Olivares C. et al., Sci. Total Environ. 599–600 (2017) 934–943;

ThP-60 / CITALOPRAM PHOTODEGRADATION: IDENTIFICATION OF TRANSFORMATION PRODUCTS BY MASS SPECTROMETRY

Rodrigo Osawa⁽¹⁾ - **Maria da Conceição Oliveira**⁽²⁾ - **Olinda Monteiro**⁽¹⁾ - **Maria Helena Florêncio**⁽¹⁾

Faculdade de Ciências, Universidade de Lisboa, Lisboa (1) - Instituto Superior Técnico, Universidade de Lisboa, Lisboa (2)

Keywords: Transformation products; time-of-flight mass spectrometry; citalopram; wastewater samples.

Introduction

Citalopram (CIT), an antidepressant drug, is an emerging contaminant frequently detected in environmental

matrices, but there are few data on its transformation products (TPs) formed in sewage treatment processes [1,2]. The objective in this study was to elucidate the formation in real matrices of CIT TPs by photodegradation under UV irradiation using LC-QTOF/MS for identification purposes.

Methods

Photodegradation experiments were conducted using a photoreactor with 100 mL of spiked distilled water (4 mg L⁻¹) cooled by water circulation. The UV irradiation source was a 450 W medium-pressure mercury-vapor lamp. The experiments consisted of 0.1 mL aliquots withdrawn in certain times. The QTOF Impact II mass spectrometer (Bruker Daltonics, Germany) was operated in the ESI positive mode in full scan and full scan product ion.

Results

The photodegradation experiments almost completely removed CIT (C₂₀H₂₁FN₂O) in 4 hours and it was possible to identify 12 probable TPs with different kinetics of formation and degradation behavior. Analyses were based on accurate mass and on the fragmentation patterns observed in the MS/MS spectra and the mass errors were less than 5 ppm. Reactions of oxidation (neutral loss of CO), hydroxylation (neutral loss of H₂O), formation of amide (neutral loss of NHCO) and dehalogenation (neutral loss of HF) were observed. The most abundant TP, namely TP4 (C₂₀H₂₂N₂O₂), was formed by the dehalogenation of fluorine and hydroxylation on the furan ring. Desmethylcitalopram (TP2, C₁₉H₁₉N₂OF) and citalopram amide (TP11, C₂₀H₂₃N₂O₂F), human metabolites, were also detected in these experiments and showed that the compounds can be formed under UV-light. A pathway for the formation of TPs was proposed. Two TPs (TP2 and TP3, C₂₀H₂₁N₂O₂F) were detected in effluent samples, where the tertiary treatment is UV disinfection.

Conclusions

LC-QTOF/MS proved to be a useful tool for elucidating TPs identify formed under laboratory conditions. UV irradiation showed good efficiency in the removal of CIT and all TPs formed completely or partially degraded until the end of the experiment (4 hours), most of them having the maximum intensity between 1 and 2 hours. Analysis of wastewater samples showed that only TP2 and TP3 were detected, presumably because in sewage treatment plant complex interactions may occur.

Novel Aspect

The results obtained are relevant to elucidate the identity of TPs formed in photodegradation and also which ones, may exhibit greater toxicity than the parent compound.

References

1. Collado N., Rodriguez-Mozaz S., Gros M., Rubirola A., Barceló D., Comas J., Rodriguez-Roda I., Buttiglieri G., Environmental Pollution, 185, 202-2012 (2014).
2. Silva L. J. G., Pereira A. M. P. T., Meisel, L. M., Lino C. M., Pena A., Science of the Total Environment, 490, 279-287 (2014).

ThP-61 / A NOVEL ANALYTICAL INTEGRATED APPROACH BASED ON MASS SPECTROMETRY FOR THE EVALUATION OF IMIDACLOPRID EXPOSURE DURING CNS DEVELOPMENT

Alice Passoni - Renzo Bagnati - Massimiliano De Paola - Enrico Davoli - Alessandro Mariani

IRCCS - Istituto di Ricerche Farmacologiche Mario Negri, Environmental Health Sciences Department, Milan

Keywords: neonicotinoid, high-resolution mass spectrometry, neurotoxic effects, multiple reaction monitoring

Introduction

The neonicotinoid pesticide imidacloprid (IMI) has been developed to act selectively as agonists of nicotinic

acetylcholine receptors (nAChRs) in insects. However, recent evidences show that it is able to bind mammalian nAChRs, leading to toxic responses in cultured neurons [1]. To properly assess the risk in humans, we developed an integrated analytical strategy combined to the analysis of biological effects in a mouse model of prenatal exposure.

Methods

The passage of blood-brain barrier (BBB) and placenta was verified analysing IMI concentrations in mouse dam brains and embryos by a multiple reaction monitoring (MRM) method. Different metabolites of IMI were identified in biological samples by high-resolution mass spectrometry (HRMS). Biological effects were evaluated in brains from newborn mice in utero-exposed to IMI by immunochemical techniques.

Results

Results of MRM analysis showed quantifiable IMI levels in both plasma and brain (about 0.3% of the total administered amount of IMI) of treated dams, suggesting a passage of mammalian BBB. HRMS analysis allowed the identification of eight main metabolites, mainly due to phase I metabolic reactions in both plasma and brain, including desnitro-IMI that has a binding affinity to nAChRs similar to nicotine. Analysis of embryos showed detectable levels of IMI (0.03%) and its metabolites, suggesting a passage through placenta. The measured concentration of IMI in biological samples has been correlated with neurodevelopmental impairments in newborn mice (verified as expression of synaptic protein and marker of immature neurofilaments in brains) and immune system alterations (reduction of microglial cells in different brain areas).

Conclusions

The newly developed analytical methods, including both MRM and HRMS, allowed an accurate quantitation of IMI in biological samples and the evaluation of relative distribution of its main metabolites. The results highlighted the passage of IMI and its main metabolites through the biological barriers (BBB and placenta) after oral exposure in a mouse model and neuro/immune system development impairments in newborn mice exposed in utero to IMI.

Novel Aspect

Combining an accurate and sensitive analytical strategy to neurological investigations at human exposure levels, our approach paves the way for a proper risk assessment.

References

1. Kimura-Kuroda J., Komuta Y., Kuroda Y., Hayashi M., Kawano H. PLoS ONE, 7, 2 (2012).

ThP-62 / NON-TARGETED APPROACH FOR THE STUDY OF LABILE DISSOLVED ORGANIC MATTER

Claudia Patriarca⁽¹⁾ - *Per J. R. Sjöberg*⁽¹⁾ - *Jonas Bergquist*⁽¹⁾ - *Lars J. Tranvik*⁽²⁾ - *Jeffrey A. Hawkes*⁽¹⁾

Uppsala University, Department of Chemistry, Uppsala (1) - Uppsala University, Department of Ecology and Genetics, Uppsala (2)

Keywords: labile dissolved organic matter, non-targeted analysis, high resolution mass spectrometry, HILIC

Introduction

Dissolved organic matter (DOM) is a complex and ubiquitous organic mixture involved in a multitude of biogeochemical processes. In order to better understand its role in the global carbon cycle, the investigation on the most elusive DOM fraction is of critical importance. The goal of this study is the development of a non-targeted method for the analysis of the labile (biologically reactive) DOM, in order to monitor its depletion and transformation on short timescales.

Methods

Labile DOM, mainly characterized by hydrophilic compounds, was produced by cold water extraction of Norway maple leaves. This material was incubated with native bacteria over four weeks. Samples were pre-concentrated and analyzed by hydrophilic interaction chromatography (HILIC) coupled to high resolution mass spectrometry (Orbitrap –MS). Data analysis and interpretation were performed by R and Matlab statistical software.

Results

The analysis of leaf-derived DOM showed a consistent transformation over the incubation period. Compounds found in the leaves exudates, characterized by high oxygen-to-carbon and hydrogen-to-carbon ratios (typically identified as sugars) were gradually removed from the solutions. Preliminary data revealed that from the transformation and depletion of the labile DOM fraction, a molecular signature typically found in the aquatic environments emerged. This suggests that recalcitrant DOM can be derived from the fast degradation of the elusive labile DOM by rapid microbial activity. Contrarily, a large part of the organic material showed no reactivity in the tested time scale suggesting that many species are not easily accessible for bacterial degradation.

Conclusions

The preliminary data show promise for the development of mechanistic studies involving the investigation of labile DOM material. The method allows the monitoring of natural compounds in the complex mixture and might be applied to find biomarker useful in environmental metabolomics analysis. Such experiments can also be used to gain insight into the global carbon cycle.

Novel Aspect

Easy and flexible non-targeted method for determination of labile DOM compounds.

ThP-63 / SCREENING OF PHARMACEUTICAL AND THEIR METABOLITES IN FISH BY HRMS USING DATA-BASE SUPPORT

Juan Manuel Peña-Herrera⁽¹⁾ - *Nicola Montemurro*⁽¹⁾ - *Sandra Pérez*⁽¹⁾ - *Damià Barceló*^(1,2)

Water and Soil Quality Research Group, Department of Environmental Chemistry (1) - IDAEA-CSIC, Barcelona, Spain (2)

Keywords: Organic pollutants, HRMS, Data-base support, Fish

Introduction

Occurrence of pharmaceuticals in the aquatic environment has raised concerns about adverse effects on exposed wildlife. Little is currently known on exposure levels of drugs in fish [1], but studies already reported detection of pharmaceuticals in fish samples [1-4]. Pharmaceuticals could be much more concentrated in fish than in polluted living waters. Additionally, fish are known to possess a hepatic detoxification system which are likely capable of metabolizing pharmaceuticals taken up from polluted waters.

Methods

Riverine fish from were collected between 2014-2017 from different regions of Spain. Bile samples were analyzed after a protein precipitation and isolation by UPLC-HRMS. Muscle were treated with organic solvents and analyzed similarly. UPLC-HRMS system present great value in establishing the fragmentation pathways to identify pharmaceuticals and elucidate the structure of fish metabolites.

Results

The HRMS data allowed screening for suspected pharmaceuticals and their metabolites and provided plausible chemical formulae. The comparison of MS/MS spectra of the parent compounds and their metabolites allowed to propose chemical structures for possible metabolites in fish bile. With this analytical methodology some metabolites, corresponding to different reactions that includes products of hydroxylation, glucuronide conjugates were identified, similar to the metabolites reported previously in our group in a similar study for carbamazepine. The suspect analysis of bile samples allowed the detection of several pharmaceuticals. Psycho-active drugs were one of the most commonly detected drugs. Their identities were proposed by matching their accurate MS and MS/MS data against different libraries.

Conclusions

The UPLC-HRMS analysis of fish bile and muscle extracts allowed the detection of pharmaceuticals and metabolite. We are presently identifying more metabolites which they still require confirmation. Psycho-active drugs have been the most commonly detected drugs in fish. Thus, this analytical methodology highlights that UPLC-HRMS is a powerful tool for qualitative analysis, allowing the search for suspected metabolites, their identification and possible quantitation of target compounds.

Novel Aspect

we propose the evaluation of the metabolism of frequently detected drugs in fish, performing a rapid screening of bile and muscle extract by HR-MS for the presence of stable intermediates.

References

1. Álvarez-Muñoz, S. Rodríguez-Mozaz, A. L. Maulvault, A. Tediosi, M. Fernández-Tejedor, F. VandenHeuvel, M. Kotterman, A. Marques, D. Barceló. *Environmental Research* 143 (2015) 56–64
2. Huerta B, Jakimska A, Gros M, Rodríguez-Mozaz S, Barceló D. *Journal of Chromatography A*. 2013;1288(0):63-72.
3. Arya, G., Tadayon, S., Sadighian, J., Jones, J., Mutsert, K., Huff T. B., & Foster, G. D. *Journal of Environmental Science and Health, Part A*, (2017) 52:7, 686-696, DOI: 10.1080/10934529.2017.1312975
4. Yi Zhang, Wen Guo, Zhenfeng Yue, Li Lin, Fengjuan Zhao, Peijin Chen, Weidong Wu, Hong Zhu, Bo Yang, Yanyun Kuang, Jiong Wang. *Journal of Chromatography B*, 1051 (2017) 41–53

ThP-64 / SUSPECTED SCREENING OF NATURAL TOXINS IN SURFACE WATER BY LIQUID CHROMATOGRAPHY HIGH RESOLUTION MASS SPECTROMETRY

Massimo Picardo^(1,2) - **Daria Filatova**^(1,2) - **Oscar Núñez**⁽²⁾ - **Marinella Farré**⁽¹⁾

Spanish Council for Scientific Research (IDAEA-CSIC) (1) - University of Barcelona, Analytical chemistry, Barcelona (2)

Keywords: Q-exactive Orbitrap, natural toxins, surface water, non-target screening, HPLC-ESI/HRMS

Introduction

Natural toxins are the prevalent contaminants in surface water reservoirs that may cause several health effects to humans [1]. Only microcystin-LR was regulated by the EU legislation [2]. High-resolution mass spectrometry (HRMS) is a powerful approach used for the non-targeted screening of natural toxins in surface water [3]. In this context, a HRMS non-targeted analysis was carried out using an in-silico approach for the identification of a wide group of biotoxins.

Methods

A mix of natural toxins representative compounds has been used as standards. After solid phase extraction (SPE) enrichment and elution, compounds were separated by a C18 HPLC column and analyzed with a ESI-

QExactive-MS. Data acquisition was based on full scan in positive ionization mode from 100 to 1500 Da with a resolution of 70000 FWHM. A homemade dataset was used as a library with Compound discoverer 2.1 software, in order to obtain a biotoxin tentative list.

Results

A double layer solid phase for SPE enrichment has been optimized for the widest range of natural toxins. The method allows good recoveries ranging from 60 to 70% and LOD below 1 µg/L. To the authors knowledge' no previous works reported an approach intended to retain such a large class of compounds. No cleanup was applied since loss of polar analytes was observed. The method was applied to perform a natural toxins suspect screening in water samples obtained from several rivers in Catalonia (NE Spain). The spectral data have been processed by Compound Discoverer 2.1 revealing the possible presence of several natural toxins in the "Riera de Rubí" river (Catalonia). Several suspected toxic compounds have been found, most of them being plant toxins and mycotoxins (but no cyanotoxins), with a level 3 of confirmation[4].

Conclusions

A novel approach able to carry out a screening of suspected natural toxins in water samples has been developed and validated. To the authors knowledge' this is the first method able to detect various groups of natural toxins simultaneously in one analysis.

Novel Aspect

This is a novel method to analyze a large group of biotoxins. A home-made dataset coupled with an in-silico approach coupled to HRMS enhanced the screening performances.

References

1. Yen, H. K. ; Lin, T. F. ; Liao, P. C. ; Toxicon, 2011, 58, 209-218.
2. Official Journal L 327 , 22/12/2000 P.0001 - 0073.
3. Schymanski, E. L. ; Singer, H. P. ; Ferguson P. L., Environmental Science & Technology, 2017, 51, 11505-11512.
4. Schymanski, E. L. ; Singer, H. P. ; Slobodnik, J. ; Ipolyi, I. M. ; Oswald, P. ; Krauss, M. ; Schulze, T. ; Haglund, P. ; Letzel, T. ; Grosse, S. ; Thomaidis, N. S. ; Bletsou, A. ; Zwiener, C. ; Ibáñez, M. ; Portolés, T. ; De Boer, R. ; Reid, M. J. ; Onghena, M. ; Kunkel, U. ; Schulz, W. ; Guillon, A. ; Noyon, N. ; Leroy, G. ; Bados, P. ; Bogialli, S. ; Stipanìev, D. ; Rostkowski, P. ; Hollender, J., Anal. Bioanal. Chem., 2015, 407, 6237.

ThP-65 / AIRBORNE MEASUREMENT OF ORGANIC GASES AND PARTICLES IN THE ATMOSPHERE WITH PTR-MS

Felix Piel^(1,2) - ***Tomas Mikoviny***⁽³⁾ - ***Markus Müller***⁽¹⁾ - ***Sven Arne Schiller***⁽²⁾ - ***Stefan Feil***⁽¹⁾ - ***Gernot Hanel***⁽¹⁾ - ***Paul Mutschlechner***⁽¹⁾ - ***Alfons Jordan***⁽¹⁾ - ***Lukas Märk***⁽¹⁾ - ***Philipp Sulzer***⁽²⁾ - ***Armin Wisthaler***⁽³⁾

Ionicon Analytik GmbH, Innsbruck (1) - University of Innsbruck, Institute for Ion Physics and Applied Physics, Innsbruck (2) - University of Oslo, Departement of Chemistry, Oslo (3)

Keywords: PTR-MS, CHARON, VOC, Organic aerosol, Airborne measurement

Introduction

Organic gases and particles are important trace constituents of the Earth's atmosphere. They affect climate and human health, and play an important role in atmospheric chemistry. Aircraft measurements of organics are analytically challenging because of the high sampling frequency requirement (> 1 Hz), low analyte concentrations in the remote atmosphere (ppt levels), and the chemical complexity of the sample.

Methods

Proton Transfer Reaction-Mass Spectrometry (PTR-MS) is an online chemical ionization mass spectrometry technique for detection of VOCs [1]. The recently developed CHARON ("CHemical Analysis of aeRosolONline") inlet enables PTR-MS instruments to measure the organic composition of submicrometer particles [2]. Significant sensitivity improvements have been recently achieved by the inclusion of radio frequency ion lenses.

Results

Our PTR-MS instruments are nowadays routinely flown on NASA's airborne science laboratories for measuring trace organics in the Earth's atmosphere. During the three most recent deployments, we i) measured a set of organic trace gases (methanol, acetone, dimethyl sulfide) that are emitted from or taken up by the North Atlantic Ocean, ii) chased another research aircraft that was fueled with a biofuel blend to measure its emissions both on the ground (during taxing) and in the air, and iii) test flew the new CHARON inlet to characterize particulate organic matter emitted from wildfires and anthropogenic sources. We will show selected data examples to illustrate how the latest generation of PTR-MS instruments has been successfully used on NASA's airborne platforms for characterizing the organic composition of gases and particles the Earth's atmosphere.

Conclusions

A PTR-MS instruments has been successfully used for airborne measurements of selected organic trace gases (*e.g.*, methanol, acetone, dimethyl sulfide) at high measurement frequency (1-10 Hz) and low detection limits (single-to-double-digit ppt at 1 second signal integration time). First test flights have been performed with the new CHARON particle inlet.

Novel Aspect

Airborne VOC measurements have been performed with unprecedented instrument sensitivity and detection limit. The CHARON particle inlet has been flown for the first time.

References

1. B. Yuan B, A. R. Koss, C. Warneke, M. Coggon, K. Sekimoto, J. A. de Gouw, Chem. Rev. 117, 13187-13229, 2017
2. M. Müller, P. Eichler, B. D'Anna, W. Tan, A. Wisthaler, Anal. Chem. 89, 10889-10897, 2017

ThP-66 / DISSIPATION STUDIES OF DIMETHACHLOR AND METABOLITES IN SOIL AND WATER BY UHPLC-ORBITRAP-MS

Rosalía López Ruíz - Roberto Romero González - José Luis Martínez Vidal - Antonia Garrido Frenich

Research Group "Analytical Chemistry of Contaminants", Department of Chemistry and Physics, Research Centre for Agricultural and Food Biotechnology (BITAL), University of Almeria, Agrifood Campus of International Excellence, ceiA3, E-04120 Almeria, Spain

Keywords: Dimethachlor, LC-Orbitrap-MS, degradation, soil, water

Introduction

Dimethachlor, which belongs to the class of chloroacetanilide, is a selective herbicide. It has a variety of metabolites, detected in soils, surface water or air [1] but they have not included in the maximum residue level of dimethachlor yet. This herbicide was not authorized in Spain, buy it was used in several countries as France, and there are scarce studies focused on the degradation pathway of this herbicide.

Methods

Extraction methods based on QuEChERS approach, using acetonitrile and magnesium sulfate and sodium

chloride as well as solid phase extraction, using Oasis HLB cartridges have been developed in order to determine dimethachlor and related metabolites in soil and water, respectively. Then, ultra high performance liquid chromatography (UHPLC) coupled to high-resolution mass spectrometry, applying Exactive-Orbitrap as analyzer, has been used.

Results

In relation to the degradation studies of dimethachlor, it has been carried out in two types of soils, sandy loam and loam, at normal and double dose and in water under sunny and darkness conditions at normal and fivefold dose. It can be concluded that it degrades slowly in both matrices. For instance, when spiked soils were monitored, the initial concentration was 35 mg/kg and in seven days, the concentration is 31 mg/kg. In relation to metabolites, a targeted analysis determines that metabolites like dimethachlor CGA, dimethachlor SYN, dimethachlor SYN Na salt and dimethachlor oxalamic acid are detected 15 days after application in soils and after 30 days in water. On the other hand, untargeted analysis can allow us to discover other metabolites like CGA 39981, CGA 369873, CGA 72649, CGA 42443 and CGA 102935.

Conclusions

A new analytical method has been developed and validated for the determination of dimethachlor and related compounds by UHPLC-Orbitrap-MS. Due to the importance of the appearance of metabolites, UHPLC-Orbitrap-MS was used in targeted mode to monitor all the compounds generated in the process and untargeted mode to look for new compounds. Degradation of the parent compound occurred slowly in soils and water.

Novel Aspect

Evaluation of the degradation of dimethachlor in soil and water applying UHPLC-Orbitrap-MS for the identification of metabolites by targeted and non-targeted approaches.

References

1. EFSA (European Food and Safety Authority). EFSA Journal. 169, 1-111 (2008).

ThP-67 / PAVING THE WAY FOR THE ANALYSIS OF TRACE COMPOUNDS IN METALLURGICAL GASES BY PROTON-TRANSFER-REACTION TIME-OF-FLIGHT MASS SPECTROMETRY. IMPORTANCE OF BACKGROUND

Jorge Iván Salazar Gómez⁽¹⁾ - *Christian Klucken*⁽¹⁾ - *Martha Sojka*⁽¹⁾ - *Holger Ruland*⁽¹⁾ - *Kevin Kähler*⁽¹⁾ - *Liudmyla Masliuk*⁽¹⁾ - *Thomas Lunkenbein*⁽²⁾ - *Robert Schlögl*^(1,2)

Max-Planck Institute for Chemical Energy Conversion, Department of Heterogeneous Reactions, Stiftstrasse 34-36, 45470 Mülheim a.d. Ruhr, Germany (1) - Fritz-Haber-Institute of the Max-Planck-Society, Department of Inorganic Chemistry, Faradayweg 4-6, 14195 Berlin, Germany (2)

Keywords: Proton-transfer-reaction time-of-flight mass spectrometry (PTR-TOF-MS), Volatile organic compounds (VOCs), coke oven gas, blast furnace gas, catalyst poisoning

Introduction

In the detailed trace analyses of industrial gases and catalytic processes the background level of carrier gases is of paramount importance. Wrong estimation may result in the underestimation of purification units [1]. Here, we focus on the analysis of the background of inert gases with different qualities at various humidity levels by proton transfer reaction (PTR) using PTR-QiTOF-MS.

Methods

Measurements were performed with a Proton-Transfer-Reaction Quadrupole interface Time-Of-Flight Mass Spectrometer (PTR-QiTOF-MS) [2]. Prior to the measurements of the background, the measurement parameters for the substances of interest were optimized, which focuses on the minimization of the

intensities of parasitic ions like N_2^+ , H^+ , NO^+ , O_2^+ and fragmentation products. By means of a commercial gas calibrator, toluene was used for the optimization process.

Results

With regard to the typical parasitic ions (N_2^+ , H^+ , NO^+ , O_2^+ , N_3^+), our results revealed that the most important effect was the change of humidity in the system. It was observed that some of these parasitic ions, such as N_2^+ and $N_2^+ H^+$ showed a significant signal reduction in the presence of humidity, whereas ions like NH_3^+ , H^+ and NO^+ remained unchanged. To date only the parasitic ions NO^+ , O_2^+ and water cluster H_3O^+ , $(H_2O)_n$ have been considered in the literature as most dominant. In our study, the most dominant parasitic ion is $N_2^+ H^+$ even after the optimization process. Our developed methodology is straightforward and allowed for assigning the origin of individual parasitic ions. Parasitic ions that are independent of the gas quality are exclusively produced in the ion source, whereas some parasitic ions which depend on humidity may be produced in the quadrupole interface. For most VOCs humidity showed an increase in sensitivity. For toluene no significant change was observed, which is in disagreement with previous results [3].

Conclusions

The characterization of VOCs in inert gases of different qualities with a PTR-QiTOF-MS and their humidity dependency in background measurements, were beneficial for the identification of instrument specific peaks. By changing the humidity it was possible to sense substances, which in dry conditions due to fragmentation were not measurable. Humidity is necessary in order to reduce the presence of parasitic ions generated in the quadrupole ion guide.

Novel Aspect

A comprehensive determination of trace contaminants present in the background of typical inert gases by means of the prototype PTR-QiTOF-MS has not been reported before. A straightforward method to identify instrument specific substances is shown.

References

1. Urban W., Lohmann H., Salazar Gómez J. I., Journal of Power Sources, 193, 359 (2009).
2. Sulzer P., Hartungen E., Hanel G., Feil S., Winkler K., Mutschlechner P., Haidacher S., Schottkowsky R., Gunsch D., Seehauser H., Striednig M., Jurschik S., Breiev K., Lanza M., Herbig J., Märk L., Märk T. D., Jordan A., International Journal of Mass Spectrometry, 368, 1 (2014).
3. Warneke C., Van der Veen C., Luxembourg S., De Gouw J. A., Kok A., International Journal of Mass Spectrometry, 207, 167 (2001).

ThP-68 / NOVEL APPLICATION OF MOLECULARLY IMPRINTED POLYMERS WITH GAS CHROMATOGRAPHY-ATMOSPHERIC PRESSURE IONIZATION-MASS SPECTROMETRY FOR HIGH THROUGH-PUT DETERMINATION OF 16 PAHS

Fereshteh Shahhoseini - Ali-Asghar Golbabanezhad-Azizi - Stefana Egli - Christina Bottaro

Memorial University of Newfoundland, Department of Chemistry, St. John's

Keywords: APGC-MS/MS, PAHs, Thin film extraction, water sample, Porous material

Introduction

Polycyclic aromatic hydrocarbons (PAHs) are common contaminants in marine and freshwaters. Identified as toxic, carcinogenic, and teratogenic, limits on the content of PAHs in water, air, and food are set by national and international regulatory agencies, e.g. the United States Environmental Agency (US-EPA). Traditional methods for extraction of PAHs from water are LLE and SPE [1,2] with excessive organic solvent use, and poor recoveries for certain species.

Methods

A thin film of molecularly imprinted polymer (MIP), designed to be selective for PAHs in water, was prepared on a glass substrate using a drop-casting method and UV-initiated radical polymerization. PAHs were extracted from 20 mL of aqueous sample into a MIP film over a fixed interval. Sorbed PAHs were extracted into a small quantity of solvent and analyzed using gas chromatography-atmospheric pressure ionization (APGC) coupled to a Xevo-TQ-S. (Waters).

Results

APGC parameter scan be optimized for soft ionization ideal for MS/MS measurements, *e.g.*, dry source conditions were used here to produce molecular ions (rather than protonated molecules), to give the highest PAH signal-to-noise. Our PAH-specific MIP films are used much like solid phase microextraction and share common approaches to method optimization. For example, sample salinity, stirring rate and extraction time, desorption solvent, and desorption time, were optimized for the 16 PAHs. The best limits of detection (LODs) were achieved when sorption to the MIP reached equilibrium at ~3 h, but since there is little sample handling simultaneous extraction of many samples allows for high throughput. Since these MIP films are intended for single-use to eliminate any carry-over and cross contamination, inter-film reproducibility was evaluated and found to be <10%. The optimized MIP-APGC-MS/MS method was applied to real samples, *e.g.*, produced water from oil and gas operations, to give fast, reproducible results at sub-ng/mL LOD.

Conclusions

A high throughput method is reported using porous thin MIP films with APGC-MS for detection of 16 PAHs from water. The methodology offers simplicity, reproducibility and speed, while minimizing consumption of organic solvents and waste generation. results showed high extraction efficiency with good LOD which meet the maximum allowable concentrations set by the USEPA.

Novel Aspect

Single use MIP thin films designed for compatibility with APGC-MS/MS are developed as adsorption material for a highthroughput sample preparation for regulated PAHs in water.

References

1. Brum, D. M. ; Cassella, R. J. ; Netto, A. D. P., *Talanta*, 74, 1392-1399 (2008)
2. Brown, J. N. ; Peake, B. M., *Analytica Chimica Acta*, 486, 159-169 (2003)

ThP-69 / NEW ANALYSIS METHOD FOR ORGANOTIN COMPOUNDS BY TANDEM-LC/MS WITHOUT THE DERIVATIZATION STEP

Jung Eun An - Yoon Chae Song - Young Soo Kim - Kyung Cheol Yi

Korea Apparel Testing and Research Institute, Environment analysis center, Anyang city

Keywords: Organotin compounds, LC/MS, Derivatization, Direct analysis

Introduction

The analysis of organotin compounds which are endocrine disrupter is important in consumer product. The derivatization step with NaBEt₄ for the analysis, however the derivatization reagent NaBEt₄ is hard to obtain because it's highly expensive, flammable and only air delivery except Europe countries. In this study, we have developed a method for direct analysis method for organotin compounds by Tandem-LC/MS without the derivatization step.

Methods

According to international environment regulatory trends, 8 kinds of organic tin compounds with one, two, or three alkyl substituents were analyzed. MBT, DBT, TBT, TPHT, MOT, DOT, TOT, and TCyT were obtained from Accustandard and Chiron. Analytical instrumentation was developed using a Tandem-LC / MS system consisting of Agilent Technologies' 1200 series HPLC and 6460 Triple Quad Mass Spectrometer.

Results

Determination of optimum extraction solvent using Design of Experiment(DoE)

Different polarity of organotin compounds should be considered when determining the best extraction solvent for simultaneous extraction of various types of organotin compounds. For example, tri-alkylated organotin with a relatively large number of alkyl groups can be regarded as a nonpolar compound compared to mono, di-alkylated organotin, and non-polar compounds such as hexane may be more suitable as the extraction solvent. Therefore, in this study, Tropolone, a complexing agent, was added to methanol used as an extraction solvent in EN ISO 17353 and ISO / TS 16179, which is a conventional organic tin compound analysis method, and used as an extraction solvent.

Determination of Tandem-LC/MS analysis condition for direct method development

Tandem-LC/MS generates ion fragments from the ion source of injected analytes using the MRM mode and transfers only the precursor ion to the collision cell in the primary mass filter (Q1) It is a system that detects only the most characteristic ions by passing through secondary mass filter (Q2) by shipping. Direct analysis of organotin compounds using Tandem-LC/MS has developed in this study does not need a derivatization process through complicated pre-treatment steps.

Conclusions

A GC/MS analysis for organotin compounds in food and consumer products are popular method, but it needs a difficult derivatization step. The new method without derivatization by LC/MS is much simpler and safe and not time consuming. The direct method is using methanol solvent containing complexing agent, tropolone. The optimum condition of extraction solvent was confirmed then the peak of each compound was obtained in chromatogram . Also international round test showed the validation of direct method by LC/MS.

Novel Aspect

Direct analysis of organotin compounds using Tandem-LC/MS does not need complicated derivatization step. Therefore the extraction step is very short, much simpler and safer operators

References

1. J. H. Banoub, J. Miller-Banoub, G. V. Sheppard, Electrospray tandem mass spectrometric measurements of organotin compounds, *Spectroscopy* 18, 95-112 (2004)
2. X. Wang, H. Jin, L. Ding, H. Zhang, H. Zhang, Organotin speciation in textile and plastics by microwave-assisted extraction HPLC-ESI-MS, *Talanta* 75, 556-563 (2007)
3. EPA Method 8323, Determination of organotins by micro-liquid chromatography-electrospray ion trap mass spectrometry

ThP-70 / CONTINUOUS-INTRODUCTION AEROSOL-TO-LIQUID SAMPLER AND ITS COUPLED USE WITH ICP-MS FOR ON-LINE SENSITIVE ANALYSIS OF ATMOSPHERIC AEROSOL

Yong-yang Su - Zhiming Li - Wei Wang - Sui Fang

Northwest Institute of Nuclear Technology, Northwest Institute of Nuclear Technology, Xi'an

Keywords: Continuous-Introduction Aerosol-to-Liquid Sampler; ICP-MS; on-line Analysis; Atmospheric Aerosol

Introduction

With the use of sampling pump [1], self-aspiration assembly [2] or aerosol introduction interface [3], ICP-MS is capable of on-line monitoring of atmospheric aerosol. However, parameters matching and conditions optimization are required prior to analysis. Moreover, concentration enrichment or matrix removal are necessary to overcome the matrix effect and sensitivity deterioration. Many types of aerosol-to-liquid samplers have been developed for efficient collection of aerosol into the liquid flow for subsequent analysis, *i.e.*, ion chromatography [4], flow cytometry [5], GC-MS [6], LC-MS [7], ESI-MS [8], novel sensor [9,10].

In order to make full use of the advantages of the aerosol sampler and ICP-MS while avoid the possible interfere, a Continuous-Introduction Aerosol-to-Liquid Sampler (CIALS) was proposed and used in couple with ICP-MS for on-line, sensitive analysis of atmospheric aerosol. This paper presents a preliminary study conducted in the laboratory.

Theory and Methods

Taking into account the arrangement of the three nozzles in the commercial BioSampler (SKC Inc.), a new collector matching the top part of BioSampler was proposed. One inlet for makeup liquid and one outlet for hydrosol are placed in the tangential but opposite position at the bottom of the collector. When CIALS is used in couple with ICP-MS, the makeup liquid flow rate (Q_1) is controlled precisely by the peristaltic pump, the hydrosol flow rate (Q_2) equals the nebulization flow rate of the ICP-MS, and the liquid loss rate (Q_3) in the collector is dependent on the sampling gas flow rate mainly due to evaporation. The sampling mechanism of this sampler is identical to the BioSampler. Because the relative position of the nozzle and the collector wall is not changed, the size-dependent collection efficient should be very close to that of BioSampler [11].

By adjusting the speed of the peristaltic pump, the liquid mass in the collector could be kept nearly constant, which is crucial for efficient sampling of aerosol into the liquid. Assuming that individual droplet is fully dissolved as soon as it's collected in the liquid, the sample concentration in the collector liquid would decay exponentially with time. The time constant equals the initial liquid mass (M_0) divided by Q_2 .

Experimental and results

The above model was verified preliminarily through experiment in the laboratory. Firstly, ICP-MS (Model 8800, Agilent Inc.) was tuned by introducing 2% HNO_3 through the sample tube and multi-elements standard solution (with concentration of around 1×10^{-9} g/g) through the "IS" tube simultaneously. The multi-elements standard solution was introduced throughout the experiment to calibrate the sensitivity of the total system. Secondly, the CIALS was connected to the nebulization system of ICP-MS through a PTFE tube. The liquid flow was continuously extracted for sequent nebulization and analysis. Meanwhile, the makeup 2% HNO_3 was introduced to compensate the liquid mass variation in the collector. The inlet of CIALS was placed downstream of a HEPA filter to sample filtered air. The background of target elements in the total system were monitored by ICP-MS when the vacuum pump was running. Thirdly, 10 microliters of standard solution containing natural Ce, Tb and Ho (214.5×10^{-6} g/g) was transferred through pipette and dropped into the collector liquid when the vacuum pump was turned off. Finally, the signals of ^{205}Tl (only from "IS" solution), ^{140}Ce (from "IS" solution and sample), ^{159}Tb and ^{165}Ho (both only from sample) were monitored by ICP-MS.

The relative variation of ^{205}Tl intensity is $\sim 7\%$, much higher than the typical value ($\sim 2\%$) when CIALS was not coupled with ICP-MS. Hence, intensity correction using internal standard is necessary. It's clearly observed that ^{140}Ce signal is a sum of stable signal from IS solution and variable signal from sample liquid. The signal profiles of ^{159}Tb and ^{165}Ho after adding of individual droplet both satisfy exponential decay law very well. The fitted values (magnitude and time constant) of the two elements math well with each other, indicating that they synchronously disperse, diffuse, and transport in the system. By using signal integration method and internal standard correction, it is possible to obtain the total amount of element in the individual particle. The fitted value of the time constant is ~ 1000 s, which roughly agrees with the value (814 s) predicted by the model. This suggests that the current system is able to monitor atmospheric aerosol with time resolution of ~ 15 min.

Conclusions and Novel Aspect

a Continuous-Introduction Aerosol-to-Liquid Sampler (CIALS) capable of on-line sampling, efficient collection of particle into liquid and continuous introduction to ICP-MS was proposed. A theoretical model for predicting the sample concentration in the collector liquid was derived. The CIALS was used in couple with ICP-MS and preliminary experiments was conducted in the laboratory. results show that the total amount of element in individual soluble particles could be determined by using the current system. The combination use of CIALS and ICP-MS has a time resolution of ~15 min for on-line monitoring of atmospheric aerosol.

References

1. M. D. Seltzer, R. B. Green, Instrumentation for continuous emissions monitoring of airborne metals, *Process Control Qual.* 6 (1994) 37–46.
2. K. Nishiguchi, K. Utani, E. Fujimori, Real-time multielement monitoring of airborne particulate matter using ICP-MS instrument equipped with gas converter apparatus, *J. Anal. At. Spectrom.* 23 (2008) 1125. doi:10.1039/b802302f.
3. Y. Y. Su, Z. M. Li, M. Li, X. W. Yi, G. Q. Zhou, J. Xu, L. H. Zhai, G. Y. Wei, F. R. Zhu, Development of Aerosol Sample introduction Interface Coupled with ICP-MS for Direct introduction and Quantitative Online Monitoring of Environmental Aerosol, *Aerosol Sci. Technol.* 48 (2014) 99–107. doi:10.1080/02786826.2013.861894.
4. D. A. Orsini, Y. Ma, A. Sullivan, B. Sierau, K. Baumann, R. J. Weber, Refinements to the particle-into-liquid sampler (PILS) for ground and airborne measurements of water soluble aerosol composition, *Atmos. Environ.* 37 (2003) 1243–1259. doi:10.1016/S1352-2310(02)01015-4.
5. D. A. Orsini, K. Rhoads, K. Mcelhoney, E. Schick, D. Koehler, O. Hogrefe, A Water Cyclone to Preserve Insoluble Aerosols in Liquid Flow — An Interface to Flow Cytometry to Detect Airborne Nucleic Acid, *Aerosol Sci. Technol.* 42 (2008) 343–356. doi:10.1080/02786820802072881.
6. J. Parshintsev, R. Räsänen, K. Hartonen, M. Kulmala, M. L. Riekkola, Analysis of organic compounds in ambient aerosols collected with the particle-into-liquid sampler, *Boreal Environ. Res.* 14 (2009) 630–640.
7. J. Parshintsev, T. Hyötyläinen, K. Hartonen, M. Kulmala, M. L. Riekkola, Solid-phase extraction of organic compounds in atmospheric aerosol particles collected with the particle-into-liquid sampler and analysis by liquid chromatography-mass spectrometry, *Talanta.* 80 (2010) 1170–1176. doi:10.1016/j.talanta.2009.09.004.
8. A. P. Bateman, S. a Nizkorodov, J. Laskin, A. Laskin, High-Resolution Electrospray Ionization Mass Spectrometry Analysis of Water-Soluble Organic Sampler †, 82 (2010) 8010–8016.
9. N. Rastogi, M. M. Oakes, J. J. Schauer, M. M. Shafer, B. J. Majestic, New Technique for Online Measurement of Water-Soluble Fe (II) in Atmospheric Aerosols, *Environ. Sci. Technol.* 43 (2009) 2425–2430.
10. F. Shen, M. Tan, Z. Wang, M. Yao, Z. Xu, Y. Wu, J. Wang, X. Guo, T. Zhu, Integrating silicon nanowire field effect transistor, microfluidics and air sampling techniques for real-time monitoring biological aerosols, *Env. Sci Technol.* 45 (2011) 7473–7480. doi:10.1021/es1043547.
11. K. Willeke, X. Lin, S. A. Grinshpun, Improved Aerosol Collection by Combined Impaction and Centrifugal Motion, *Aerosol Sci. Technol.* 28 (1998) 439–456. doi:10.1080/02786829808965536.

ThP-71 / ACETONE DEGRADATION IN A PLASMA REACTOR STUDIED IN REAL-TIME BY PTR-MS

Sebastien Thomas⁽¹⁾ - **Michel Heninger**⁽¹⁾ - **Nicole Blin-simiand**⁽²⁾ - **Lionel Magne**⁽²⁾ - **Hélène Mestdagh**⁽¹⁾ - **Stéphane Pasquier**⁽²⁾ - **Essyllt Louarn**⁽¹⁾

Laboratoire de chimie physique, Univ. Paris Sud, CNRS, Univ. Paris Saclay, 91405, Orsay (1) - Laboratoire de physique des gaz et plasmas, Univ. Paris Sud, CNRS, Univ. Paris Saclay, 91405, Orsay (2)

Keywords: CIMS, volatile organic compound, FTICR, PTRMS, pollutant removal

Introduction

Acetone is one of the major pollutants in the Earth's atmosphere [1]. So that it is becoming more and more necessary to reduce acetone anthropogenic emission. The new non-thermal plasma technique for pollution control has received increasing attention [2]. While there has been recent progress in acetone removal, clear identification of the reaction products, eventually toxic, to be done.

Methods

We used a photo-triggered low-temperature plasma discharge [3] to study the mechanisms of acetone decomposition at different concentrations in N_2 and N_2/O_2 mixtures. Degradation by-products are detected in-situ by a compact FTICR mass spectrometer BTrap [4], that uses chemical ionisation by proton transfer reaction from H_3O^+ [5] or by charge transfer from O_2^+ . Identification and direct quantification of the molecules is conducted in real-time in the instrument.

Results

Presence of oxygen in the gas mixture reduces the discharge efficiency, even for a small O_2 amount, up to 5%. The decomposition of the molecule entails the formation of many by-products. For acetone diluted in N_2 , at least 20 different compounds were detected. The major detected by-products include nitrile compounds (hydrocyanic acid and acetonitrile), formaldehyde and acetaldehyde, alcohols (mainly methanol and ethanol). Some products, such as hydrogen cyanide, appear at the first discharge and increase rapidly. These products are then reprocessed but not completely. These results are compared with those of acetone decomposition in N_2/O_2 , which leads to oxygenated compounds such as nitro-alkanes.

Conclusions

The large number of by-products leads to better understanding of the mechanisms involved in plasma degradation. Our work shows that the non-thermal plasma technique degrades efficiently acetone in N_2 and N_2/O_2 gas mixtures. We plan to extend these analysis methods to others discharge conditions with or without associating catalysts.

Novel Aspect

Coupling of a high resolution compact FTICR to a VOC degradation plasma reactor and real-time monitoring of by-products.

References

1. H. Singh, D. Ohara, D. Herlth, J. Geophys. Res., 99, 1805-1819 (1994)
2. A. Vanderbrouke, R. Morent, N. De Geyter, C. Leys, J. Hazard. Mat., 195, 30-54 (2011)
3. B. Lacour, V. Puech, S. Pasquiers, Recent Res. Dev. Appl. Phys., 6, 149-191 (2003)
4. H. Latappy, J. Lemaire, M. Heninger, E. Louarn, E. Bauchard, H. Mestdagh, Int. J. Mass Spectrom., 405, 13-23 (2016)
5. A. Ellis, C. Mayhew, Proton Transfer Reaction Mass Spectrometry: Principles and Applications (2014), John Wiley & Sons: Chichester, UK.

ThP-72 / UPTAKE OF DIOCTYL PHTHALATE BY POPULUS ALBA VILLAFRANCA

Francesca Vannucchi⁽¹⁾ - ***Erika Carla Pierattini***⁽¹⁾ - ***Andrea Raffaelli***⁽²⁾ - ***Alessandra Francini***⁽¹⁾ - ***Luca Sebastiani***⁽¹⁾

Biolabs, Institute of Life Science, Scuola Superiore Sant'Anna, Pisa (1) - Istituto di Fisiologia Clinica, CNR, Pisa (2)

Keywords: DOP, micro-pollutants, phthalates, poplar, xenobiotics

Introduction

Phthalates are micro-pollutants detected in aqueous systems and their removal through wastewater treatment plants results ineffective [1,2]. The accumulation and metabolism of organic pollutants (caffeine, erythromycin, sodium dodecyl sulfate) in poplar has been demonstrated [3-6] and in this study we evaluated the poplar capability to uptake and translocate dioctyl phthalate inside the plant.

Methods

Plants of *Populus alba* Villafranca clone, growing in a hydroponics condition, were treated for 21 days with 0, 40 and 400 $\mu\text{g L}^{-1}$ of dioctyl phthalate (d_4 -DOP). The health status of the plants was monitored through the plant growth, relative chlorophyll content and fluorescence analyses. The d_4 -DOP concentrations in roots, stem and leaves were evaluated through LC-MS/MS analysis.

Results

Treatments did not affect significantly poplar health at none concentrations tested: plants maintained a good photosynthetic performance and no variation amongst relative chlorophyll contents were observed between control and treated d_4 -DOP plants. After 21 of treatments, 400 $\mu\text{g L}^{-1}$ d_4 -DOP induced an increase in root dry biomass. On the contrary, leaves number and leaf dry weight were reduced. Poplar revealed the ability to uptake the d_4 -DOP with the highest accumulation observed in roots. Indeed, the d_4 -DOP concentration in roots reached $7.1 \pm 3.29 \mu\text{g g}^{-1}$ FW and 3.4 ± 2.28 , after 1 day and 21 days of treatments respectively. In leaves d_4 -DOP concentration values were lower than $0.02 \mu\text{g g}^{-1}$ FW.

Conclusions

P. alba Villafranca showed a good tolerance to d_4 -DOP at both concentrations tested and the accumulation of d_4 -DOP occurred mainly in roots after one day of treatment.

Novel Aspect

For the first time tree removal of DOP from aqueous solutions was explored. Our findings can contribute to the development of green technologies for the remediation of phthalates.

References

1. Abdel Daiem M. M., Rivera-Utrilla, J., Ocampo-Pérez, R., Méndez-Díaz, J. D., Sánchez-Polo, M. *Journal of environmental management*, 109, 164-178 (2012).
2. Luo Y., Guo W., Ngo H. H., Nghiem L. D., Hai F. I., Zhang J., Liang S., Wang X. C. *Science of the Total Environment*, 473, 619-641 (2014).
3. Pierattini E. C., Francini A., Raffaelli A., Sebastiani L. *Environmental Science and Pollution Research*, 23, 7298-7307 (2016).
4. Pierattini E. C., Francini A., Raffaelli A., Sebastiani L. *Science of the total Environment*, 569–570, 540–547 (2016).
5. Pierattini E. C., Francini A., Raffaelli A., Sebastiani L. *Tree physiology*, 38, 109-118 (2018).
6. Raffaelli A., Pierattini E. C., Francini A., Sebastiani L. *Comprehensive Analytical Chemistry*, 241-266 (2018).

ThP-73 / IDENTIFICATION OF DEGRADATION PRODUCTS OF SELECTED PHARMACEUTICALS IN SOILS BY HIGH RESOLUTION MASS SPECTROMETRY

Andrea Vojs Staňová^(1,2) - **Katerina Grabicova**⁽¹⁾ - **Adam Borik**⁽¹⁾ - **Radka Kodesova**⁽³⁾ - **Ales Klement**⁽³⁾ - **Roman Grabic**⁽¹⁾

University of South Bohemia in České Budějovice, Faculty of Fisheries and Protection of Waters, South Bohemian Research Center of Aquaculture and Biodiversity of Hydrocenoses, Vodňany, Czech Republic (1) -

Department of Analytical Chemistry, Faculty of Natural Sciences, Comenius University in Bratislava, Bratislava, Slovak Republic (2) - Czech University of Life Sciences Prague, Faculty of Agrobiological Sciences, Department of Soil Science and Soil Protection, Prague 6, Czech Republic (3)

Keywords: mass spectrometry, pharmaceuticals, degradation products, soils

Introduction

Pharmaceuticals are frequently detected in different matrices, including soils, in the last decade [1, 2]. Different degradation/transformation products can be formed by interaction of soil bacteria with pharmaceuticals [3]. Only limited information on this topic is available, we aimed to study the fate of the selected pharmaceuticals and their possible metabolites in the different soils.

Methods

Fortified soils samples were incubated under controlled conditions, extracted and analysed with LC-ESI-full scan HRMS (QExactive, Thermo Fisher Scientific) at a resolution of 70,000 FWHM. Compound Discoverer 2.0 software was used for data processing. The high resolution product scan (HRPS) at 17500 FWHM resolution was applied to list of suspect metabolites to obtain the MS2 fragmentation.

Results

The degradation of 6 pharmaceuticals (carbamazepine, citalopram, clindamycin, fexofenadine, irbesartan and sulfamethoxazole) was investigated in 4 different soil types. The parent compounds degraded in all of the soils. Dissipation curves for the parent compounds of interest were constructed to evaluate a degradation rate by means of slope values. We obtained list of suspect transformation products from two different workflows in CD software for each compounds. The list was further evaluated using HRPS data to confirm or disprove proposed compounds. Tentatively identified TPs were semiquantified with method described by Koba [4]. Time trends of TPs concentration could further help in both degradation mechanism and kinetics clarification

Conclusions

Using sequence of full scan HRMS, advanced data processing, HRPS structure confirmation and back semi quantification from HRMS we were able to identify degradation products of selected pharmaceuticals and obtained important information on transformation of pharmaceuticals in soil.

Novel Aspect

Identification and semi quantification of degradation products of pharmaceuticals in soils significantly increase information value of analysis.

References

1. Kodesova R., Grabic R., Kocarek M., Klement A., Golovko O., Fer M., Nikodem A., Jaksik O., Sci. Total Environ., 511, 435-443 (2015).
2. Kodešová R., Kočárek M., Klement A., Golovko O., Koba O., Fér M., Nikodem A., Vondráčková L., Jakšík O., Grabic R., Sci. Total Environ., 544, 369-381(2016).
3. Koba O., Golovko O., Kodešová R., Fér M., Grabic R., Environmental Pollution 220, 1251-1263 (2017).
4. Koba O., Steinbach C., Kroupová H. K., Grabicová K., Randák T., Grabic R., Rapid Commun. Mass Spectrom., 30, 1153-1162 (2016)

Acknowledgements

This study was supported by the by the Grant Agency of Czech Republic (No. 17-8937S)

ThP-74 / THE STUDY OF SAMPLING PIPELINE EFFECT FOR GASES PARTIAL PRESSURE MEASURE WITH MASS SPECTROMETER

Rongxin Yan

Beijing Institute of Spacecraft Environment Engineering, CAST, Beijing, China

Keywords: Mass spectrometer, Gases partial pressure measurement, Sampling pipeline, Stability time

Introduction

In some gases analysis or partial pressure measure with mass spectrometer, The gases are far from the mass spectrometer because mass spectrometer is not portable or the gases environment is toxic. The sampling pipeline connecting gases and mass spectrometer has to be used. But the sampling pipe how to affect measure results such as stable time, ingredient aberration, which are worthwhile to study.

Methods

Based on the study of unsteady-state gas flow with the effect of pipeline adsorption, the physical mechanism of gases at molecular regime in pipeline flowing into mass spectrometer chamber and pumped out by vacuum system has been studied. The gases transportation equation in sample pipeline and gases partial pressure distributing equation in mass spectrometer chamber are set up. The effect of sample pipeline for some gases partial pressure in air are calculated.

Results

The relation equation of sample pipeline gases pressure, the mass spectrometer chamber gases pressure, analysis time, pumping rate to the chamber, gas molecular mass, pipeline length has been got.

The pressures in sampling pipeline gas are a complex exponential equation with sampling time and sampling pipeline distance;

The gases balance pressures in mass spectrometer chamber are direct proportion to the gases pressure of the sampling place;

The sampling pipeline change the ration of gases ingredient. After sampling equilibrium, the ratio of partial pressure is inverse to the square root of the ratio of gas mass. The pressure in mass spectrometer chamber is inverse to the length of pipeline.

The minimum pressure stable time in chamber is proportion to the square of pipeline length and inverse to gas molecule thermal motion speeds.

Novel Aspect

The mass spectrometer gas pressure equation with sample pipeline length and sample time was deduced. The stable time and ingredient aberration of gases partial pressure measure with mass spectrometer can be calculated.

ThP-75 / DIRECT ANALYSIS OF 17 PERFLUORINATED COMPOUNDS IN WATER AT LOW PARTS-PER-TRILLION LEVELS BY LC-MS/MS WORKFLOW

Sha Josh Ye⁽¹⁾ - *Jingcun Wu*⁽¹⁾ - *Feng Qing*⁽¹⁾ - *Frank Kero*⁽²⁾ - *Luca Piatti*⁽²⁾ - *Stefan Edler*⁽²⁾ - *Derek J. Mattern*⁽²⁾ - *Marco Gori*⁽²⁾

Perkin Elmer Canada, Discovery & Analytical Solutions, Woodbridge (1) - Perkin Elmer, Discovery & Analytical Solutions, Waltham (2)

Keywords: Direct Analysis, LC-MS/MS, perfluorinated compounds

Introduction

Perfluorinated compounds (PFCs) are used as surfactants, fire-retardants, nonstick cookware coatings, and coatings for paper packaging for over a half century[1,2]. PFCs have recently received attention as emerging persistent organic pollutants [3-14]. Presently, a LC-MS/MS direct sample injection method was developed for analysis of PFCs in water samples. The results indicate that a good sensitivity, robustness were achieved

with LOQ below the proposal by EU.

Methods

The PFCs native and internal standards were obtained from Wellington Laboratories (Guelph, Ontario). 100 μ L of lake or tap water samples were first filtered then injected onto a PerkinElmer Brownlee SPP C18 (100x2.1 mm, 2.7 μ m) column for separation at 0.5 mL/min with total run time of 15 minutes. The mobile phases were A) water and B) methanol, both containing 5 mM ammonium formate. Detection was via a QSight triple quadrupole LC-MS/MS system.

Results

The calibration standards were prepared with different levels of PFCs to determine instrument detection limits and set up calibration curves. Detection limits for most analytes were ≤ 1 ng/L except for PFBA, PFTA, PFHxDA, and PFODA which were ≤ 5 ng/L. The calibration curves for all analytes were linear over three orders of magnitude with a R^2 value of > 0.995 . Identification and confirmation of the analytes in investigated samples was achieved by comparing their respective retention times and peak area ratios between the quantifier and qualifier MRM transitions with those from the standards. Ten (10) samples were screened for PFOA and PFOS with the developed method. Five of them were detected with PFOA at a concentration of larger than 2 ng/L and six of the samples were detected with PFOS at a concentration of larger than 5 ng/L.

Future work would include further development of shorter LC separation time, robustness evaluation, and the analysis of other brands of bottle water and water from rivers or lakes using the workflow developed in this study.

Conclusions

The fast and simple LC-MS/MS method for 17 perfluorinated compounds analysis in water samples was developed by coupling a UHPLC system to a QSight 220 triple quadrupole mass spectrometer. This method can be employed to the direct determination of PFCs in drinking and surface water samples, with LOQs well below the limits set by regulatory boards from many countries.

Novel Aspect

Direct analysis of perfluorinated compounds in water at low ppt levels with improved sensitivity and robustness by a simple LC-MS/MS workflow.

References

1. Schultz, M. M. ; Barofsky, D. F. ; Field, J. A. Environ. Eng. Sci. 20, 487-501 (2003).
2. U. S. Environmental Protection Agency Public Docket AR226-1486. : Washington, DC (2003).
3. Boulanger, B. ; Peck, A. M. ; Schnoor, J. L. ; Hornbuckle, K. C. Environ. Sci. Technol. 39, 74-79 (2005).
4. Giesy, J. P. ; Kannan, K. Environ. Sci. Technol. 35, 1339-1342 (2001).
5. Kannan, K. ; Corsolini, S. ; Falandysz, J. ; Fillmann, G. ; Kumar, K. S. ; Loganathan, B. G. ; Mohd, M. A. ; Olivero, J. ; Van Wouwe, N. ; Yang, J. H. ; Aldous, K. M. Environ. Sci. Technol. 38, 4489-4495 (2004).
6. Schultz, M. M. ; Barofsky, D. F. ; Field, J. A. Environmental Engineering Science, 20, 487-501 (2003).
7. Ellis, D. A. ; Mabury, S. A. ; Martin, J. W. ; Muir, D. C. G. Nature 412, 321-324 (2001).
8. PUBLIC HEALTH STATEMENT Perfluoroalkyls, Public Health Service Agency for Toxic Substances and Disease Registry (2015)
9. PERFLUORINATED CHEMICALS (PFCs): PFOA & PFOS INFORMATION PAPER, Association of State and Territorial Solid Waste Management Officials (2015)
10. ACT SHEET PFOA & PFOS Drinking Water Health Advisories, US Environmental Protection Agency (2016).
11. News on limit values proposed for PFASs, metals. Chemical Watch (2018).
12. McCleaf P, Englund S, Östlund A, Lindegren K, Wiberg K, Ahrens L. Water Research. 120, 77-87 (2017).
13. Villagrasa, M. ; de Alda, M. L. ; Barcelo, D. Anal. Bioanal. Chem. 386, 953-972 (2006).
14. van Leeuwen, S. P. J. ; de Boer, J. Journal of Chromatography A 1153, 172-185 (2007).

ThP-76 / MOLECULAR COMPOSITION OF DISSOLVED ORGANIC MATTER IN THE TALDICE ICE CORE

Roberta Zangrado ⁽¹⁾ - ***Veronica Zanella*** ⁽²⁾ - ***Ornela Karroca*** ⁽²⁾ - ***Elena Barbaro*** ⁽¹⁾ - ***Natalie M. Kehrwald*** ⁽³⁾
- ***Dario Battistel*** ⁽²⁾ - ***Andrea Gambaro*** ⁽²⁾ - ***Carlo Barbante*** ⁽¹⁾

Consiglio Nazionale delle Ricerche, Institute for the Dynamics of Environmental Processes, Mestre-Venezia (1) - Ca' Foscari University of Venice, Department of Environmental Sciences, Informatics and Statistics, Mestre-Venezia (2) - U. S. Geological Survey, Geosciences and Environmental Change Science Center, Denver (3)

Keywords: Untargeted analysis, ice cores

Introduction

Deep polar ice cores are archives that collect information on paleoclimatic and paleoenvironmental conditions during glacial-interglacial swings. This information is then obtained by detailed analyses of the ice to obtain past temperature, greenhouse gas concentrations, and sea ice extent etc. The presence of trace chemicals in ice, reflect changes in their sources, transport, deposition, preservation.

Methods

An untargeted analysis of dissolved organic matter (DOM) provides a general overview of molecular species in the ice samples giving information on their sources and processes. For this kind of analysis typically 100 - 900 mL are required, this is not possible when dealing with deep ice cores. Using a nano-UPLC-nano-ESI-HRMS technique we have detected the major molecular species using only 4 μ L of melted ice.

Results

This work allowed us to tentatively identify generally aliphatic species, that were attributable mainly to the presence of saturated and unsaturated fatty acids, hydroxyl fatty acids and their degradation products, along with species from the oxidation of isoprene and monoterpenes and other compounds more difficult to interpret. Interglacial had more species while their number decreased in glacial samples. The reduction in oxidation products of isoprenes and monoterpenes is probably due to temperature changes that affected terrestrial vegetation and the sea ice extent. Passing from interglacial to glacial times, unsaturated FAs almost disappeared because they are exposed for longer to oxidative processes.

Conclusions

Until now untargeted analyses have been limited to surface snow and shallow ice cores covering a short time span. To the best of our knowledge these are the first indications of the behavior of fatty acids and their oxidation products and the oxidation products of isoprene and monoterpenes in a deep ice core.

Novel Aspect

The innovative aspect of this work was the creation of a micro-volume analytical method that allowed the molecular characterization of DOM in a deep ice core.

ThP-77 / INVESTIGATIONS ON THE ADSORPTION OF REES BY GO WITH DIFFERENT SIZES: MECHANISM AND APPLICATION

Wenbo Zhao

Chinese Academy of Geological Sciences, National Research Center for Geoanalysis, Beijing

Keywords: Rare earth elements (REEs); (Inductively coupled plasma mass spectrometry) ICP-MS; Graphene oxide (GO); Adsorption

Introduction

Rare earth elements (REEs) are major resources for modern society but have the risk of environmental pollution. With a high surface area, carbon nanomaterial shows the good adsorption character and has the bright future to become a novel REEs separation material [1]. Herein, we demonstrate the special adsorption behavior on REEs by graphene oxide. The mechanism is also investigated and has the chance of being applied in the comprehensive ecological improvement.

Methods

GO nanosheets were synthesized from natural flaky or expanded graphites using the modified Hummer methods. The physicochemical properties of GO nanosheets were characterized by SEM, FTIR and Raman respectively. We tested the adsorption behavior of GO nanosheets on Rare earth elements with different concentrations. After centrifugation and filtration, the solution was transferred into ICP-MS to obtain residual content of REEs.

Results

We obtained GO with 5 different meshes from natural flaky graphites and 3 giant GO with pre-expanding process. The SEM results exhibited that GO surface morphology was well nanoscale thin layer. FTIR spectra indicated that GO had similar O-containing functional groups as mentioned [2]. It was worth mentioning that GO of 500 and 200 meshes had higher carboxyl, aromatic, sharp C=O and –OH, where as the giant GO of 100 mesh had higher –OH. Raman spectroscopy showed that GO of 300 and 100 meshes had more defects/edges. Through the key ICP-MS detection, we found GO had superior adsorption ratio on REEs to other inorganic elements even including heavy metals. Then, we compared the adsorption behavior on 17 REEs with 5 different concentrations. When at proper concentration of REE, GO with 500,200 and 100 meshes exhibited more excellent adsorption ratio. Surprisingly, neither a lower nor a higher concentration of REEs could display adsorption difference. The accordance of adsorption and characterization made us to illuminate the adsorption mechanism properly.

Conclusions

Comparing with other inorganic elements, we verified the superior adsorption on REEs by GO which was fabricated by modified Hummers methods. The corresponding adsorption mechanism was illuminated that aromatic, carboxyl, sharp C=O and –OH contribute to the whole special adsorption behavior. The proper GO meshes should have to be selected to separate or enrich REEs in the future environmental pollution.

Novel Aspect

GO has special adsorption behavior on REEs superior to other inorganic elements and the mechanism is also been illuminated.

References

1. Wang T., Zhao D., Guo X., Correa J., Riehl B., Heineman W., Analytical Chemistry, 86, 4354-4361(2014).
2. Chen X., Chen B., Environmental Science & Technology, 49, 6181-6189(2015).

ThP-78 / HILIC-MS (NON-)TARGET SCREENING, MONITORING AND TOXICOLOGICAL EVALUATION OF HIGHLY POLAR COMPOUNDS IN DRINKING WATER AND SOURCES

Annemieke Kolkman⁽¹⁾ - ***Rosa Sjerps***⁽²⁾ - ***Pascal Kooij***⁽¹⁾ - ***Kirsten Baken***⁽²⁾ - ***Pim de Voogt***⁽²⁾

KWR Watercycle Research Institute, Laboratory for Materials Research and Chemical Analysis, Nieuwegein, Netherlands (1) - KWR Watercycle Research Institute, Chemical Water Quality and Health, Nieuwegein, Netherlands (2)

Keywords: very polar compounds, non-target screening, hydrophilic interaction liquid chromatography, Orbitrap Fusion

Introduction

Highly polar organic compounds are of special interest for drinking water utilities, since these substances are likely to end up in drinking water [1]. Currently there is an analytical gap, a monitoring gap, and a lack of toxicity data for persistent and mobile organic compounds [2]. We aimed to close these gaps by applying a target and non-target screening using hydrophilic interaction liquid chromatography (HILIC) coupled to high resolution mass spectrometry in a monitoring study [3], followed by a toxicological risk assessment.

Methods

HILIC-MS target and non-target screening approach including a pre-treatment was developed. The method was validated for the determination of 38 highly polar compounds in surface water and in drinking water. 29 samples from surface water, river bank filtrate, ground water and drinking water in The Netherlands and Flanders were analysed.

A toxicological risk assessment was performed for the observed highly polar compounds.

Results

More than half of the 38 highly polar target compounds were present in surface-, ground and/or drinking water at levels between LOD and LOQ or above the LOQ. Traces of several highly polar compounds were observed in drinking water. One compound, dichloroacetic acid, was introduced during disinfection treatment.

Non-target HILIC-MS screening detected 145 features: compounds with a unique combination of accurate mass and retention time. The identity of eleven compounds could be confirmed by Reference standards, the remaining compounds are yet to be identified. Due to threshold settings for detection not all target compounds that were observed with target HILIC analysis were observed with non-target HILIC-MS screening.

A pGLV could be derived for 8 chemicals. For the remaining compounds present in the evaluated samples, insufficient information was available to derive a provisional drinking water guideline value (pGLV).

Conclusions

A HILIC-MS non-target screening approach was developed and validated for the determination of highly polar compounds. The HILIC target analysis and non-target screening were proven useful for the detection of highly polar compounds in source waters and in drinking water. Further identification of chemicals observed and toxicological risk assessment are imperative to evaluate the relevance of highly polar chemicals in (sources of) drinking water.

Novel Aspect

An novel non-target screening approach using HILIC was developed for highly polar compounds revealed the presence of known and unknown polar compounds in the water cycle.

ThP-79 / DIRECT ANALYSIS OF PHENOLIC COMPOUNDS IN CORK WASTEWATER BY HPLC-QTOF

Francisco Javier Yuste Córdoba - Belén Godoy Cancho

CICYTEX, ICMC, Mérida, Spain

Keywords: cork, phenolic compounds, HPLC, QTOF

Introduction

Cork is the outer surface of the cork oak tree and it is a versatile material that can be used for a variety of products, mainly stoppers. Industrial preparation of cork requires the immersion of cork planks for approximately one hour in boiling water. Cork boiling wastewaters (CBWs) are waste streams with high loads of phenolic compounds (1). This work studies the direct analysis of CBWs with HPLC/QTOF in order to

identify useful by-products for their recovery.

Methods

The samples were analysed by HPLC on a reversed-phase column RP-18 (250 nm x 4 mm x 5 µm). and detection was carried out using a diode array detector. Other sub-samples were analysed by HPLC-QTOF in order to confirm the target compounds.

Results

Several low-weight phenols were detected with HPLC-DAD by comparing UV-vis absorption maxima in mass spectral analysis and the retention times of standards. Mass spectra of these phenolic compounds were compared with the corresponding standards by QTOF. 18 target compounds were studied: gallic, protocatechuic, ellagic, p-hydroxyphenyl acetic, p-hydroxybenzoic, vanillic and syringic acids; esculetin, vanillin, syringaldehyde, ferulic, caffeic, p-coumaric and salicylic acids; coniferyl aldehyde, sinapaldehyde, eriodictyol and naringenin. From these phenolic compounds, significative concentrations of vanillin and gallic, ellagic, vanillic, protocatechuic, ferulic and syringic acids were detected, confirmed and quantified. Besides, some characteristic ellagitannins were also analyzed but they were not detected.

Conclusions

The direct analysis of CBWS by HPLC/QTOF confirmed several low weight phenolic compounds highlighting the presence of vanillin and gallic, ellagic, vanillic, protocatechuic, ferulic and syringic acids. This kind of compounds and their concentrations coincides in general terms with other studies (1). Regarding extractable compounds from cook such as ellagitannins didn't appear, maybe due to the higher temperature of the boiling process.

Novel Aspect

This work shows that the direct analysis of phenolic compounds in CBWs without pre-concentration achieved good results and made the tests simpler.

References

1. I. P. Marques, L. Gil, F. La Cara, *Biotechnology for Biofuels*, 7, 67 (2014).

ThPS-S05 - Safe & Valuable Food

Chairs: Michel W.F. Nielen, Emanuela Gregori

ThP-80 / A NEW ANALYTICAL METHOD FOR THE QUANTIFICATION OF THREE LIGNANS IN ESPRESSO COFFEE BY USING HPLC-MS/MS TRIPLE QUADRUPOLE

Simone Angeloni⁽¹⁾ - ***Giovanni Caprioli***⁽¹⁾ - ***Gulzhan Khamitova***⁽¹⁾ - ***Luciano Navarini***⁽²⁾ - ***Gianni Sagratini***⁽¹⁾
- ***Sauro Vittori***⁽¹⁾

University of Camerino, School of Pharmacy, Camerino (1) - illycaffè S. p. A., Research & Scientific Coordination, Trieste (2)

Keywords: Espresso Coffee, EC, Lignans, HPLC-MS/MS

Introduction

Lignans, phenolic compounds that are present in a number of foodstuffs, possess some biological activities *i.e.*, antioxidant and antiestrogenic [1,2]. The aim of our work was to develop a new analytical method to

quantify three lignans, *i.e.*, secoisolariciresinol (SECO), lariciresinol (LARI), matairesinol (MAT) from espresso coffee (EC) by using HPLC-MS/MS Triple Quadrupole.

Methods

Different sample preparation methods, including dilutions, acidic hydrolysis and enzymatic digestions were evaluated and, once validated, the most performing was applied for the quantification of target compounds in coffee samples. The HPLC-MS/MS Triple Quadrupole, equipped with electrospray ionization (ESI) source, operating in negative ionization mode, was used for the analysis. Detection was performed in multiple reaction monitoring (MRM) mode.

Results

Each lignan showed different transitions with different intensities of MS/MS acquisition parameters in MRM mode; the highest product ion abundances for all three compounds were found at 350°C as temperature of the drying gas of ionization source. At this temperature the precursor ions for lignans were the deprotonated molecules $[M-H]^-$ in negative polarity; therefore, those conditions were chosen for lignan monitoring. About the sample preparation process, the most performing was the enzymatic digestion with the Clara-Diastase (mixture of enzymes), because it showed the highest recovery values (93.55-97.62 %). The analytical method was characterized by short chromatography run time (the separation of target compounds was obtained within 4 minutes) and high sensitivity (LODE ranging from 2 to 3 $\mu\text{g L}^{-1}$ and LOW ranging from 5 to 10 $\mu\text{g L}^{-1}$).

SACO and LAIR were found in all nine EC samples from 27.9 to 52.0 $\mu\text{g L}^{-1}$ and from 5.3 to 27.8 $\mu\text{g L}^{-1}$, respectively, contrary to MAT that it was not possible to detect it in each type of coffee.

Conclusions

A new analytical method for the quantification of three lignans that avoided long clean up and enrichment sample procedures has been developed and validated. The application of our method to EC samples showed that SECO was the lignan present at the highest concentration (27.9-52.0 $\mu\text{g L}^{-1}$) followed by LARI (5.3-27.8 $\mu\text{g L}^{-1}$). Therefore, regular consumption of EC can contribute to the dietary intake of lignans.

Novel Aspect

The developing of a new analytical method for the quantification of three lignans in espresso coffee by using HPLC-MS/MS and their determination in different coffee samples.

References

1. Landete J. M., Food Research International, 46, 410 (2012).
2. Adlercreutz H., Mazur W., Annals of medicine, 29, 95 (1997).

ThP-81 / CATABOLIC MARKERS 24,25-DIHYDROXYVITAMIN D2/D3 IN BOVINE AND PORCINE SERA SHOW INTER-SPECIES DIFFERENCES BETWEEN VITAMIN D2 AND VITAMIN D3 METABOLISM

Kirsten G. Dowling⁽¹⁾ - **George L. J. Hull**⁽¹⁾ - **Aoife Hayes**⁽¹⁾ - **Sarah K. Duffy**⁽²⁾ - **Alan K. Kelly**⁽²⁾ - **John V. O'Doherty**⁽²⁾ - **Kevin D. Cashman**⁽¹⁾

Cork Centre for Vitamin D and Nutrition Research, School of Food and Nutritional Sciences, University College Cork, Cork (1) - School of Agriculture and Food Science, University College Dublin, Dublin (2)

Keywords: 24,25-dihydroxyvitamin D, vitamin D biofortification, beef, pork, LC-MS/MS

Introduction

Supplementation of animal feeds with vitamin D2 (VitD2) is less effective than vitamin D3 (VitD3) at increasing serum 25-hydroxy vitamin D [25(OH)D; status marker] of beef cattle and pigs [1,2]. A faster rate

of degradation of VitD2 may explain this finding, and measurement of the serum 24,25-dihydroxy catabolic markers, 24,25(OH)2D2 and 24,25(OH)2D3, by LC-MS/MS would test this hypothesis.

Methods

Sera were obtained from animal feeding trials in which beef cattle received diets containing 0 IU/kg VitD or 4000 IU/kg as VitD2 or VitD3, and pigs received diets containing 2000 IU/kg as VitD2 or VitD3 (n=8/10 per group). [1,2] Sera were analysed with a validated LC-MS/MS method for 24,25(OH)2D3, 25(OH)D3 and 25(OH)D2 in 150 µL serum [3] which was extended to include the analyte 24,25(OH)2D2.

Results

The majority of the 30 beef sera analysed had concentrations of 24,25(OH)2D2 (90%) or 24,25(OH)2D3 (77%) below the LOQs (0.68 and 1.09 nmol/L, respectively). By contrast, 7 out of 8 pigs receiving VitD2 had quantifiable 24,25(OH)2D2, ranging from 1.85-4.77 nmol/L, and 24,25(OH)2D3 concentrations in pigs from both treatment groups ranged from 2.6-24.2 nmol/L. These differential 24,25(OH)2D findings were despite the fact that mean+SD serum 25(OH)D was higher in beef cattle compared to pigs (122+37 and 53+18 nmol/L, respectively), which suggests a slower 24-hydroxylation rate in beef cattle. In relation to catabolism of the vitamins amongst pigs, the mean+SD serum 24,25(OH)2D3 of pork receiving VitD3 was 21.0+2.8 nmol/L with no detected 24,25(OH)2D2, and the mean+SD serum 24,25(OH)2D3 and 24,25(OH)2D2 of pork receiving VitD2 was 4.5+1.8 nmol/L and 2.3+1.3 nmol/L, respectively.

Conclusions

Differences in serum 24,25(OH)2D2 and 24,25(OH)2D3 concentrations of beef cattle and pigs shows that regulation of the vitamin D degradation pathway is not the same in these two species. Amongst pigs, the catabolic rate of 25(OH)D2 is lower than that of 25(OH)D3, suggesting that different rates of catabolism do not account for the difference in efficacy of VitD2 and VitD3 feed supplementation.

Novel Aspect

Most 24,25-dihydroxyvitamin D methods measure only 24,25(OH)2D3. Including 24,25(OH)2D2 allowed direct comparison of vitamin D2 and D3 metabolism.

References

1. Duffy, S. K., O'Doherty, J. V., Rajauria, G., Clarke, L. C., et al., Food Chem, 256, 18-24 (2018).
2. Duffy, S. K., Kelly, A., Rajauria, G., Jakobsen, J., et al., Meat Sci (in press).
3. Dowling, K. G., Hull, G., Sundvall, J., Lamberg-Allardt, C., et al, J Chrom A, 1497, 102-109 (2017).

ThP-82 / EVALUATION OF OXIDATIVE STRESS MARKERS IN HUMAN BIOLOGICAL MATRICES BY MEANS OF UHPLC-MS/MS

Federico Fanti⁽¹⁾ - **Flaminia Vincenti**⁽²⁾ - **Nicola Bernabò**⁽¹⁾ - **Mauro Serafini**⁽¹⁾ - **Dario Compagnone**⁽¹⁾ - **Manuel Sergi**⁽¹⁾

University of Teramo, Faculty of Bioscience and Technology for Food, Agriculture and Environment, Teramo
(1) - Sapienza University of Rome, Department of Chemistry, Rome (2)

Keywords: Oxydative Stress Markers, Oxysterols, Isoprostanes, UHPLC-MS/MS

Introduction

Oxidative stress (OS) is a risk factors for human health, and several stress markers have been proposed for the evaluation of OS, and many of them are available for this purpose. The induction of organ function disorders can lead to age-related decline in biological functions in human [1]. The aim of work is to create suitable methods for identification of oxidative stress markers by means of UHPLC-MS/MS techniques

applied in human samples, in particular for Oxysterols and Isoprostanes [2-3].

Methods

In this project, Methods for extraction from various tissue has optimized for little amounts of samples (for example sperm and plasma) or for analyte enrichment (urine). After extraction, a clean-up procedure has performed in the samples using SPE methods for each type of analytes. For each of them a an UHPLC-MS/MS experiments was performed, being careful to take every possible parameter into account, using different source (ESI or APCI) and optimizing every instrument parameter.

Results

For the analytes different experiments was performed, for the oxysterols the best condition for analyst was obtained through using of APCI source operating in positive ion scan, in according with literature [4], using C18 column for the separation of analytes. In other hand, isoprostanes were more responsive by using of ESI source operating in positive ion mode; for the liquid chromatographic separation a C18 column was used. SPE clean up methods were optimized for the selected classes of analytes: for isoprostanes SPE step was used to concentrate the OS markers present in urine before injection in the auto sampler. The extraction protocols of Oxysterols has been optimized, starting from bovine sperm procedures [5]: the protocols were miniaturized to allow the use of small sample amount. The developed methods were tested for different biological samples (human).

Conclusions

UHPLC-ESI(APCI)-MS/MS methods for the detection and quantitation of OS markers were optimized. These markers were detected and quantified in different human biological matrices both for oxysterols and isoprostanes; both groups of samples were taken from subjects that followed controlled feedings, adapted to modify the levels of oxidative stress.

Novel Aspect

New analytical methods for detection and quantitation of OS markers were developed and validated, increasing sensitivity and reliability in the different matrices.

References

1. Yoshikawa T., Naito Y., Kondo M., Neuroscience, 16, 603-612 (1990).
2. Iuliano L., Micheletta F., Natoli S., Ginanni Corradini S. et al., Analytical Biochemistry, 312, 217–223 (2003).
3. Montuschi P., Barnes P. J., and Roberts II L. J. et al., The FASEB Journal, 18, 1791-1800 (2004).
4. Jiang ., Sidhu R., Porter F. D., Yanjanin N. M. et al., Journal of Lipid Research, 52, 1435-1445 (2011)
5. Brouwers J. F., Boerke A., Silva P. F. N., Garcia-Gil N. et al., Biology of Reproduction 85, 128–136 (2011)

ThP-83 / COMPREHENSIVE CHARACTERIZATION OF MINIMAS, A NATURAL COMPLEX FOOD SUPPLEMENT, BY DIFFERENT MASS SPECTROMETRY TECHNIQUES

Enrico Flamini⁽¹⁾ - *Giada Fodaroni*⁽¹⁾ - *Sara Tamimi*⁽¹⁾ - *Stella Bedont*⁽¹⁾ - *Denise Decarli*⁽¹⁾ - *Michela Burico*⁽¹⁾ - *Anna Gaetano*⁽¹⁾ - *Luisa Mattoli*⁽¹⁾ - *Pietro Traldi*⁽²⁾

Aboca Spa Società Agricola, Area Ricerca Fitochimica, Sansepolcro (1) - Città della Speranza, Istituto di Ricerca Pediatrica, Padova (2)

Keywords: MiniMas, mass spectrometry, food supplement

Introduction

Food supplements are intended to support, maintain or optimize physiological parameters, therefore the organism's homeostasis. The complex mechanisms that contribute to the maintenance of homeostasis have brought researchers to move from a reductionist to a holistic vision of physiology, biology and

medicine. In this panorama, natural complex food supplements play a key role as all the compounds presents forms a natural complex system ready to interact at various level with the complexity of human organism[1].

Methods

To reach the global characterization of Minimas, different chromatographic methods based on mass spectrometry have been used. Volatile organic compounds have been analysed by means GC-QqQ and head space sampling. Organic polar compounds have been analysed by means UHPLC-qToF. The metal part of inorganic compounds have been analysed by means ICP-MS.

Results

The quality assurance of all natural complex food supplements is a very important goal particularly that of Minimas, a food supplement that helps to maintain the homeostasis of blood pressure and cardiac functions. The old fashion of analyze only few marker compounds is not enough, as the activity of Minimas, as well that of all natural complex products, cannot be referred to the presence of few compounds. All the compounds of the product are important as they express an effective complex natural system. Therefore, a comprehensive product characterization is necessary to assure batch compliance and efficacy through the time. In this work, plant dried-extracts ingredients of Minimas and Minimas itself have been analyzed by means of GC-QqQ and UHPLC-qToF methods. Inorganics have been characterized by means of ICP-MS. All the data collected have been matched with an in-house library containing more than 1000 natural compounds providing the identification of several compounds, that than have been quantified [2].

Conclusions

The potency of the developed analytical approach applied to the natural complex food supplement Minimas led to its comprehensive characterization, getting a quali-quantitative identity card of its complexity. The analysis performed by different mass spectrometry techniques permits to check the presence of all Minimas ingredients, opening new perspective to the routine quality control of natural complex food supplements.

Novel Aspect

A complete mass spectrometry based analytical platform is available and it represents a very important tool to get the broader compositional knowledge of all natural complex products, included Minimas.

References

1. Fang F. C., Casadevall A., Infection and Immunity, 79,1401, (2011)
2. Wu H., Guo J., Chen S., Liu X., Zhou Y., Zhang X., Xu X., J Pharm Biomed Anal, 91, 267 (2013).

ThP-84 / THE SYNERGISM OF MALDI LIFT-TOF/TOF MS AND ESI-Q-TOF MS INSTRUMENTS TO DISCOVER NEW ISOFORMS OF ALLERGENS FROM HAZELNUT OLEOSIN EXTRACT

Sterfano Nebbia⁽¹⁾ - ***Cristina Lamberti***⁽¹⁾ - ***Marcello Manfredi***⁽²⁾ - ***Emilio Marengo***⁽³⁾ - ***Laura Cavallarini***⁽¹⁾ - ***Maria Gabriella Giuffrida***⁽⁴⁾

CNR, ISPA, Grugliasco (TO) (1) - ISALIT, Spin-off of DISIT, University of Piemonte Orientale, Italy, Alessandria (2) - University of Piemonte Orientale, Department of Sciences and Technological Innovation, Alessandria (3) - CNR, ISPA, Colletterto Giacosa (TO) (4)

Keywords: MALDI-TOF/TOF, ESI-Q-TOF, hazelnut allergens, allergens

Introduction

Oleosins are proteins associated to oil bodies of seeds and investigated because demonstrated to be

responsible of severe allergic reaction in allergic patients [1]. During the characterization of hazelnut oleosins, selected by immunoreaction with the sera of allergic Italian children, the analyses with two MS instruments improved our capacity of data interpretation and allowed to identified a new isoform probably more reactive than the two oleosins already known.

Methods

A protocol for the extraction of oil bodies associated proteins was optimized and SDS-PAGE was performed. The PMF method ("in gel" trypsin digestion, peptide collection and MS analysis) was applied to the bands at 14 and 17 kDa. All protein bands were identified by both MALDI-TOF/TOF and LC-MS/MS analyses. MALDI-TOF/TOF spectra were manually interpreted, if necessary, and results were used to refine LC-MS/MS interpretation.

Results

The electrophoresis profile of the oil bodies associated proteins showed two most abundant bands. MALDI TOF/TOF analysis was performed on bands at 14 and 17 kDa from SDS-PAGE of hazelnut extract, where 100% of protein content resulted to be oleosins Cora 12 and 13, already known as allergens. The LC-MS/MS analysis confirmed the same results. However, the manual checking of major signals of the MALDI-TOF spectra of the band at 17kDa revealed that only minor signals identified the oleosin Cora 12 and that at least two high intensity signals were without any protein attribution. The fragmentation spectra from LIFT-TOF/TOF of two major picks were analyzed using MS-Tag Protein Prospector software in identity mode and didn't give any protein identification. When the software was set in homology mode with Cora 12 as template, a new N-terminal was revealed indicating the presence a new isoform not present in the NCBI database.

Conclusions

In protein identification the dependence of MS techniques from the DNA and protein database is very strong and new isoforms are rarely identified. In this work a new isoform of oleosins from hazelnut was identified combining manual and software-assisted interpretation of MALDI LIFT-TOF/TOF spectra.

Novel Aspect

Sometimes old practices like to have look at spectra and try to assign almost all signals by manual interpretation of MS/MS spectra allow to really discover something new that normally miss out on automatic data interpretation.

References

1. Zuidmeer-Jongejan L, Rivas FM, Winter MG, Akkerdaas JH, Summers C, Lebens A, Knulst AC, Schilte P, Briza P. *Clinical and Translational Allergy*, 4, 4 (2014).

ThP-85 / DETERMINATION OF PYRROLIZIDINE ALKALOIDS IN PLANT MATERIAL USING ON-LINE SPE COUPLED TO UHPLC-MS/MS

Anja Grüning - Julia Sander - Sven Vedder

Shimadzu Europa, CIPS, Duisburg

Keywords: Pyrrolizidine alkaloids, on-line SPE, UHPLC-MS/MS

Introduction

Pyrrolizidine alkaloids (PAs) are secondary plant metabolites that are supposed to be carcinogenic and genotoxic. They occur mainly in plants of the Boraginaceae, Asteraceae and Fabaceae families. Plant food and beverage, phytopharmaceuticals or even animal feed can easily be contaminated with PAs and enter the food chain.

Currently there are discussions on possible regulatory measures caused by the presence of PAs in honey, tea, herbal infusions and food supplements.

Existing methods include laborious sample preparation, *e.g.* solid-liquid extraction followed by solid phase extraction for clean-up.

Here we report an on-line SPE UHPLC-MS/MS method, which overcomes the difficulties of combining

Methods

Tea samples were extracted twice with 0.05 M sulfuric acid by sonication. Before centrifugation the pH of the extracts was adjusted with ammonium hydroxide. These samples were put into the autosampler and transferred to the on-line SPE column using an aqueous solution. After washing the sample was eluted and trapped into a loop. By switching the loop the eluted sample was transferred to the analytical column. A binary gradient separated the PAs for quantification. Due to this hardware set-up UHPLC with high backpressure and on-line SPE which is pressure limited were successfully combined. By careful fine-tuning of the SPE elution and the chromatographic conditions the separation of critical peak pairs could be maintained.

Results

By using the reported instrument set-up, analysis and thus the quantification of 16 PAs (Indicine, Trichodesmine, Senkirkine, Senecivernine, Senecyphilline, Senecionine, Retrosine, Monocrotaline, Lycopsamine, Lasiocarpine, Jacobine, Intermedine, Heliotrine, Europine, Erucifoline, Echimidine) and 14 of their related *N*-Oxides (Echimidine-*N*-oxide, Lasiocarpine *N*-oxide, Jacobine *N*-oxide, Intermedine *N*-oxide, Heliotrine *N*-oxide, Erucifoline-*N*-Oxide, Europine *N*-Oxide, Indicine *N*-oxide, Senecivernine *N*-oxide, Senecyphilline *N*-oxide, Senecionine *N*-oxide, Retrosine *N*-oxide, Monocrotaline-*N*-oxide, Lycopsamine *N*-oxide) could be performed. Calibration curves showed good precision and accuracy and even in a complex matrix like tea we were able to easily quantify the PAs in at least the range of 10 to 300 µg/kg. This is consistent to established methods using manual sample preparation, but less time consuming and less tedious. For all analytes, linear weighted regression resulting in r^2 0.99 could be achieved, with $S/N > 10$ for LLOQ levels.

Conclusions

A method for high-sensitivity analysis was successfully developed for PA analysis in plant material. The manual sample preparation could be reduced to a minimum as the set up of on-line SPE followed by UHPLC-MS/MS saves additional clean-up steps without compromising the performance of the assay.

Novel Aspect

On-line SPE of pyrrolizidine alkaloids coupled to an UHPLC-MS/MS system which overcomes the problems with the different pressure ranges.

ThP-86 / FREE FATTY ACIDS AND HEXANAL AS SHELF LIFE MARKERS IN CHEESE: EVALUATION OF FOOD PACKAGING EFFICIENCY THROUGH A QUANTITATIVE HS-SPME-GC-MS METHOD

Astride Franks Kamgang Nzekoue - Elisa Vittori - Sauro Vittori - Giovanni Caprioli - Gianni Sagratini

University of Camerino, School of Pharmacy, Camerino

Keywords: Cheese shelf-life, free fatty acids, hexanal, food packaging, HS-SPME-GC-MS

Introduction

This study aims to develop a quantitative HS-SPME-GC-MS method to monitor concomitantly the evolution of short and medium chain free fatty acids and hexanal during cheese storage. Indeed, food storage is associated with the production of these markers whose amounts determine the taste and the spoilage level

of cheese [1,2]. Therefore, their analysis can be useful in the evaluation of dairy product shelf-life [3].

Methods

Hexanal, butanoic, isovaleric, hexanoic, octanoic and decanoic acids were analyzed by optimizing the GC-MS conditions in SIM mode. Moreover, different salting-out agents and SPME-fibers were compared to improve their extractions. This optimized method was validated and applied on two cheese samples (Grana Padano and Fresh Asiago) to compare the efficacy of two packaging systems: a poly coupled packaging named Ideabrill® and a HDPE film.

Results

The trapping abilities of 4 fiber types (100 µm PDMS, 75 µm CAR/PDMS, 50/30 µm DVB/CAR/PDMS and 85 µm PA) was evaluated and followed by the comparison of 4 salts (NaCl, Na₂SO₄, NaH₂PO₄ and Na₂CO₃) on the extraction of FFAs and hexanal.

The 75 µm CAR/PDMS fiber and the use of NaH₂PO₄ as salting-out agent resulted to be the most appropriate conditions for the extraction of analytes. A progressive increase of hexanal and FFAs during storage was noted confirming that these compounds are sensitive markers of cheese deterioration. All along storage, the concentrations of these markers became higher in the HDPE film. After 12 weeks of study for Grana Padano (GP) and 4 weeks for Fresh Asiago (FA), the Ideabrill® packaging allowed a lower increment of FFAs (11.6-fold increase for GP and 5.9 for FA) in comparison with the HDPE film (21.5-fold increase for GP and 7.4 for FA). The amounts of hexanal were also lower in the Ideabrill® packaging (0.58 against 1.04 mg/kg in GP and 0.38 against 0.55 mg/kg in FA).

Conclusions

By slowing down the increment of FFAs and hexanal, the Ideabrill® packaging has proved to be more effective than a common HDPE film in the preservation of the cheeses studied.

This method can represent an important solution in the food industry sector, whose growth requires the optimization of fast and performing analytical methods that could help for example in shelf life studies, quality control or active packaging development.

Novel Aspect

Hexanal and FFAs are good chemicals markers of cheese shelf-life. This is the first GC-MS method reported to quantify simultaneously these markers in cheese.

References

1. McSweeney, P. L. H., Sousa, M. J., Lait, 80, 293-324 (2000).
2. Mannion, D. T., Furey, A., Kilcawley, K. N., International Journal of Dairy Technology, 69, 1-12 (2016).
3. Elisia, I., Kitts, D. D., Journal of Clinical Biochemistry and Nutrition, 49, 147-152 (2011).

ThP-87 / BIOACTIVE COMPOUNDS OF BRAZILIAN GREEN COFFEE BEANS COLLECTED IN DIFFERENT LEVELS OF MATURATION BY ESI-LTQ-ORBITRAP

Mayara Fumiere Lemos^(1,2) - ***Consuelo Perez***⁽²⁾ - ***Demian Ifa***⁽²⁾ - ***Aymbire Fonseca***⁽³⁾ - ***Rodrigo Scherer***⁽¹⁾

Pharmaceutical Sciences Graduate Program, University of Vila Velha, Espírito Santo, Brazil (1) - York University, Centre for Research in Mass Spectrometry, Toronto, Canada (2) - The Brazilian Agricultural Research Corporation, Incaper (3)

Keywords: Conilon coffee; chlorogenic acids; caffeine; trigonelline

Introduction

Harvesting coffee beans at different cultivating conditions may influence their chemical composition and

the quality of the beverage [1,2]. Considering the economic importance of conilon coffee, the objective of the study was to identify the presence of metabolites in inedited conilon coffee clones and to compare them with known species, using a high resolution technique, as well as to verify if different maturation times can affect the composition of the coffee.

Methods

Green coffee beans were collected in Espirito Santo, Brazil (3 cultivars of Arabica coffee, 3 clones of conilon coffee and one cultivar of Robusta Tropical). The samples were harvested in 3 different levels of ripeness (60%, 80% and 100% of mature) and the extracts were prepared with methanol. The extracts were injected into a mass spectrometer (LTQ-Orbitrap Elite) equipped with an electrospray ionization (ESI) source, operating in positive and negative ion modes.

Results

A total of 21 coffee bean extracts were characterized in positive and negative ion modes by ESI-MS and tandem MS. Coffee bean metabolites, such as caffeine, trigonelline, and some sugars were observed in positive and a wide range of polyphenols such as chlorogenic acids in negative ion mode. It was found that coffee bean species Arabica, Robusta and clones 101,105 and 108 showed similar metabolomic profiles but distinct relative abundances. In positive ion mode, trigonelline was the most abundant compound with higher concentrations in Arabica samples and in clones 105 and 108. No trend was observed regarding to the maturation of this compound. However, the content of caffeine was higher in Robusta coffee and also in clones 105 and 108, increasing with maturation. In negative ion mode, caffeoylquinic acid was the most abundant compound with higher concentrations in Robusta coffee and in clones 105 and 108. The polyphenols and chlorogenic acids detected showed constant concentrations during the maturation in almost the samples.

Conclusions

Among the 3 clones studied, it was possible to conclude that clone 105 and 108 were similar to Arabica coffee (better quality) and Robusta (higher productivity), being promising for a good quality. Trigonelline and caffeoylquinic acid were the main compounds detected in coffee beans. It was observed that maturation may affect some metabolites and it is important to monitor these changes because these compounds are scent and flavor precursors of the final beverage.

Novel Aspect

To understand the conilon coffee quality, as well as the effects of different harvesting conditions, which could impact the composition and quality of the beverage.

References

1. Arruda N., Hovell A., Rezende C., Freitas S., Couri S., Bizzo H., Química Nova, 34, 819-824 (2011).
2. Lemos M., Lemos M., Pacheco H., Endringer D., Scherer R., Industrial Crops and Products, 70, 41-47 (2015).

ThP-88 / DEVELOPMENT OF A LC-MS/MS METHOD FOR THE DETECTION OF SPECIES-SPECIFIC MUSCLE PEPTIDES IN PAPS

Daniela Marchis⁽¹⁾ - **Marta Leporati**⁽¹⁾ - **Alessandra Altomare**⁽²⁾ - **Giancarlo Aldini**⁽²⁾ - **Cristina Casalone**⁽³⁾

Istituto Zooprofilattico Sperimentale del Piemonte, Liguria e Valle d'Aosta, S. S. Igiene dei Mangimi, Torino (1) - Università degli Studi di Milano, Dipartimento di Scienze Farmaceutiche, Milano (2) - Istituto Zooprofilattico Sperimentale del Piemonte, Liguria e Valle d'Aosta, S. C. Neuroscienze con annesso CEA, Torino (3)

Keywords: PAPS, peptides, proteins, LC-MS/MS, BSE

Introduction

Ruminant PAPs (Processed Animals Proteins) were forbidden in feed since 2001 [1], while non ruminant PAP were reauthorized in aquaculture in 2013 [2]. A PCR method is used in official controls: it permits to discriminate between ruminant and non-ruminant materials [3]. Mass-spectrometry can be helpful in identifying peptides from specie-specific muscular proteins, since milk is an allowed ingredient in feed and PCR is not able to discriminate the origin of DNA.

Methods

High resolution mass spectrometry was applied to identify specie-specific peptides from bovine muscular proteins. Then, a LC-MS/MS method was developed in order to detect these peptides in animal feed samples eventually contaminated/adulterated with bovine PAPs. Sample preparation consists of a proteins extraction step, followed by an accurate extract purification by means of SDS-PAGE and by an enzymatic proteins digestion in order to obtain the peptides [4].

Results

Three peptides, TR-12, TR-14 and YK-16, have been identified as markers of bovine PAPs presence by means of high resolution mass spectrometry. They derived from three bovine skeletal muscle proteins, Desmin, Vimentin and Myoglobin respectively. Proteins' specificity has been carefully evaluated, excluding all proteins also contained in milk, colostrum and in other sources. Then, an accurate selection has been executed on peptides found, excluding all peptides not specie-specific and with no adequate physical/chemicals properties.

A LC-MS/MS method for the simultaneous detection of these peptides in animal feed samples was developed. This method is able to detect a contamination of animal feed with bovine PAP 0.1% w/w.

Validation as qualitative method is in progress. The method was also applied to a real sample resulted non-compliant for ruminant DNA at PCR analysis and not reporting milk or milk by-products as ingredients: peptide TR-12 has been correctly identified.

Conclusions

The bovine peptides identified appear to be excellent biomarkers for detection of bovine PAP in animal feed. Their presence is an indisputable signature of bovine PAP, as they were selected from bovine meat protein sequences, and no possible cross reactivity appear to exist with dairy products. This could be of primary importance, as it will allow official and in house laboratories to detect unquestionably the presence of bovine PAP in feed.

Novel Aspect

A new approach to detect bovine PAP contamination was developed. It permits to discriminate between milk and bovine PAP addition in feed, unlike currently official methods.

References

1. European Commission, Official Journal of European Union, L147, 1-40 (2001).
2. European Commission, Official Journal of European Union, L21, 3-16 (2013).
3. European Commission, Official Journal of European Union, L20, 33-43 (2013).
4. Marchis D., Altomare A., Gili M., Ostorero F., Khadjavi A., Corona C., Ru G., Cappelletti B., Gianelli S., Amadeo F., Rumio C., Carini M., Aldini G., Casalone C., Journal of Agricultural and Food Chemistry, 65, 10638-10650 (2017)

ThP-89 / METABOLIC SURVEY OF ITALIAN VEGETABLES: FTICR-MS AND NMR-BASED UNTARGETED ANALYSES

Alessandro Maccelli⁽¹⁾ - **Anatoly P. Sobolev**⁽²⁾ - **Donatella Capitani**⁽²⁾ - **Barbara Chiavarino**⁽¹⁾ - **Simonetta Fornarini**⁽¹⁾ - **Luisa Mannina**⁽¹⁾ - **Maria Elisa Crestoni**⁽¹⁾

Sapienza Università di Roma, Dipartimento di Chimica e Tecnologie del Farmaco, Roma, Italy (1) - CNR, Istituto di Metodologie Chimiche, Laboratorio di Risonanza Magnetica "Annalaura Segre", Monterotondo, Italy (2)

Keywords: ESI FT-ICR, NMR, metabolomics, White Celery, Red Pepper

Introduction

Untargeted analyses of complex biological mixtures such as vegetable extracts is a rapidly expanding field. Combined use of Fourier transform ion cyclotron resonance mass spectrometry (FT-ICR MS) and nuclear magnetic resonance (NMR) spectroscopy allows to increase metabolome coverage devoid of authentic standards[1] and is applied here to characterize white Celery of Sperlonga (WCS) and "Cornetto di Pontecorvo" red pepper (CPRP) samples.

Methods

Both vegetables are endowed with Protected Designation of Origin and Protected Geographical Indication mark. Extracts from fruit (peel, pulp and seeds of CPRP) and plant (heart, stalk and leaves of WCS) were obtained by the Bligh-Dyer method. High resolution mass analyses were conducted with an ESI source coupled to a 4.7 T FT-ICR (BioApex, Bruker Daltonics) mass spectrometer; NMR spectra were recorded on a Bruker Avance 600 spectrometer equipped with a multinuclear z-gradient 5 mm probe head.

Results

Use of both NMR spectroscopy and FT-ICR MS has allowed a comprehensive metabolite profile of pulp, peel and seeds hydroalcoholic and organic extracts of CPRP grown in open field (OF) and in green house (GH), revealing the presence of sugars, amino acids, organic acids, polyphenols, flavonoids. The very low content of biogenic amines and mycotoxins verified quality and freshness of the vegetables. Interestingly, the main features of CPRP are preserved in either OF or GH cultivation methods. [2] The same experimental approach applied to hydroalcoholic extracts of WCS revealed that leaves and stalks present a similar composition, including choline, sugars, polyalcohols, amino acids, phenols, glycosylated compounds, free organic acids, including chlorogenic and ascorbic acids. In organic extracts, major peaks correspond to phtalides, as senkyunolide, responsible for organoleptic properties, chlorophyll a, linoleic and linolenic fatty acids, the latter more abundant in leaves.

Conclusions

The high sensitivity, resolution and mass accuracy of FT-ICR MS and the powerful quantitative and structural characterization allowed by NMR spectroscopy have provided a detailed metabolic description of different vegetable extracts. The large variety of secondary metabolites reflects the expression of biosynthetic pathways influenced by the environment.

Novel Aspect

These data may support efforts in preserving biodiversity and in developing new cultivation strategies.

References

1. Pan Z., Raftery D., Anal Bioanal Chem, 387, 525-527 (2007)
2. Sobolev A. P. et al. Food Chem, 255, 120-131 (2018).

ThP-90 / GC/MS MONITORING OF VOLATILES IN GRAPPA DURING AGING IN TWO DIFFERENT WOOD BARRELS

Mirko de Rosso⁽¹⁾ - Christine M. Mayr⁽²⁾ - Roberto Castagner⁽³⁾ - Roberto Sonogo⁽³⁾ - Riccardo Flamini⁽¹⁾

Council for Agricultural Research and Economics, Viticulture & Oenology (CREA-VE), Conegliano (TV) (1) - University of Padova, Department of Agronomy, Food, Natural Resources, Animals and Environment (DAFNE), Legnaro (PD) (2) - Distillerie Acquavite S. p. A., Visnà di Vazzola (TV) (3)

Keywords: Grappa, GC-MS, volatiles, barrel aging

Introduction

The name "Grappa" is reserved exclusively for the grape pomace distillate obtained from raw materials produced in Italy and distilled and bottled in plants located in Italy (CE Regulation n. 110/2008). This research project aimed to investigate different conditions (wood type and alcohol percentage during aging) in the production of two high quality distillates.

Methods

Grappa produced from a white marc (Prosecco) and a red blend (Cabernet Sauvignon/Merlot), at 55% and 68% Ethanol each, were aged over 6 month in 225 L, middle toasted cherry and oak barrels, which were selected on the basis of previous findings [1,2]. After one, three, and six month, samples were taken, volatiles isolated by SPE and analyzed by GC/MS. In particular, the evolution of variety aroma compounds from the marc and volatiles released from the woods, was monitored.

Results

The analysis of the different grappa samples before ageing allowed the identification of more than 50 volatiles, mainly higher alcohols, aldehydes, acids, ethyl esters and acetates. The main differences between the two types of samples was, as expected, the higher abundance of variety aroma compounds, mainly terpenes, in the white marc grappa. Samples analyzed after one month showed already a slight increase of some phenolic compounds from the wood, but due to the low levels the impact on the overall aroma can be neglected. After three month the aroma of the grappa was already characterized by a wood aged aroma, and the concentration of phenolic compounds increased significantly.

Conclusions

At the end of the process the different distillates showed characteristic profiles due to the different varieties and the wood type used. Main volatiles detected in the samples markedly characterized the organoleptic profile of the final products.

Novel Aspect

This is a systematic GC/MS investigation of the impact of wood type and alcohol percentage on the aroma profile of grappa aimed to make new products.

References

1. De Rosso M., Cancian D., Panighel A. Dalla Vedova A., Flamini R, Wood Sci. Technol., 43, 375-385 (2009).
2. Flamini R., Dalla Vedova A., Cancian D., Panighel A., De Rosso M, J. Mass Spectrom., 42, 641-646 (2007).

ThP-91 / SENSITIVE AND SPECIFIC ALLERGEN SCREENING ANALYSIS USING LC-MS/MS

Ashley Sage⁽¹⁾ - Jianru Stahl-Zeng⁽²⁾ - Chris Borton⁽³⁾ - Leesun New⁽³⁾

SCIEX, Lakeside Drive, Warrington, UK (1) - SCIEX, Darmstadt Germany (2) - SCIEX, Redwood City, CA, USA (3)

Allergens in food can result in severe or fatal reactions. At present, there is no known cure for an allergenic reaction and the only thing a person can do is to avoid the potential cause *e.g.* nuts, milk etc. Food allergens are of increasing interest due to food allergy recalls that have doubled in recent years. To help safeguard consumers from food related allergies, warning labels on packaging and in restaurants are a must

these days to allow people to make a decision on what they eat. However, food testing is also important and as a result, it is vital to have a robust and specific analytical method to reliably identify and quantitate allergens that may be present. This work presents data from a method that has been highly characterized and verified to determine several different types of allergens using a tryptic digest and LC-MS/MS analysis to measure allergen peptides with a high degree of flexibility, specificity but also with high sensitivity. Samples were tryptically digested and the resulting extracts analyzed using LC-MS/MS. Peptides that were identified to be associated with a specific allergen were measured using electrospray ionization and scheduled MRM data acquisition. In this presentation a comparison between traditional methods and the LC-MS/MS will be described. Important considerations for method development will be discussed. The ability of the method to identify allergens and to quantify them in food samples will be presented and discussed. The method results from several different food matrices will be presented to demonstrate the potential of this method.

Novel Aspect

This work demonstrates the process for identification of the correct peptide for allergen identification.

ThP-92 / LC-MS/MS: A USEFUL TOOL TO INVESTIGATE MARINE LIPOPHILIC TOXINS IN MUSSELS

Melania Siracusa - Simone Bacchiocchi - Sara Evangelista - Arianna Piersanti

Centro Referenza Nazionale controllo microbiologico e chimico Molluschi bivalvi, Istituto Zooprofilattico Sperimentale Umbria e Marche "Togo Rosati", Perugia

Keywords: Marine Lipophilic Toxins, LC-MS/MS, mussels, phytoplankton, profiles

Introduction

Marine Lipophilic Toxins (MLTs) are compounds produced by marine microalgae. They often accumulate in mussels representing a risk for human health. MLTs include okadaic acid OA, dinophysistoxins DTXs, pectenotoxins PTXs, yessotoxins YTXs and azaspiracids AZAs and are regulated within Europe [1, 2]. The Istituto Zooprofilattico Sperimentale Umbria e Marche "Togo Rosati" is responsible for their monitoring in mussels harvested in the Marche Coast (Italy).

Methods

The official LC-MS/MS method [3] was adopted to evaluate MLTs contamination in mussels bred and collected along the coast, biweekly or monthly, since March 2012 until December 2017. Mass experiments were performed using a hybrid triple-quadrupole/linear ion trap 3200 Q TRAP spectrometer (AB Sciex, Darmstadt, Germany) in Multiple Reaction Monitoring (MRM) mode. Two transitions were selected for each molecule to allow reliable quantification and identification.

Results

The data collected about MLTs in shellfish samples (*M. Gallo provincialis*) showed different toxin profiles over the years. OA and YTXs were the only groups present in the analysed mussels. The concentrations of both of them usually increase sharply in autumn-winter, starting from the Northern areas of the coast and afterwards extending southern. Until 2013 YTXs were predominant reaching concentrations that exceed the regulatory limit (MRL) of 1 mg YTX equivalents Kg⁻¹, subsequently increased to 3.75 mg YTX equivalents Kg⁻¹ [4]. In 2015 and 2017, OA was the most abundant toxin measured in the samples with levels four/five times the regulatory limit (160 µg OA equivalents Kg⁻¹), causing closures of local mussel farms. The different toxins profiles characterizing mussels in the years are very likely a result of variable climatic and environmental conditions such as temperature, rainfall, salinity, that may influence the blooms of different phytoplanktonic species.

Conclusions

The LC-MS/MS method is suitable to monitor MLTs contamination in shellfish. The analytical technique used allows to investigate the profiles in terms of abundance and variety. The hypothetical correlation between climatic and environmental parameters and the observed levels could be useful in the prediction of mussels contamination, therefore a efficient tool to evaluate this threat for public health.

Novel Aspect

Mass spectrometry applied to MLTs analysis provided detailed information on mussels contamination, enabling the study of temporal trends and toxin profiles.

References

1. Regulation (EC) N° 853/2004, Official J. Eur. Union, L139, 55–205 (2004).
2. Regulation (EC) N° 2074/2005, Official J. Eur. Union, L338, 27-59 (2005).
3. EU-RL-MB- AECOSAN, EU-Harmonised Standard Operating Procedure for determination of Lipophilic marine biotoxins in molluscs by LC-MS/MS, vers. 5 (2015).
4. Bacchiocchi S., Siracusa M., Ruzzi A., Gorbi S., Ercolessi M., Cosentino M. A., Ammazalorso P., Orletti R., *Toxicon*, 108, 115-125 (2015).

This work was funded by the Italian Health Ministry - Project Code: RC0072016

ThP-93 / ADVANCES IN MASS SPECTROMETRY FOR FOOD SAFETY AND AUTHENTICITY

Despina Tsipi⁽¹⁾ - ***Helen Botitsi***⁽²⁾

General Chemical State Laboratory, Food Safety & Quality Laboratories, Athens (1) - General Chemical State Laboratory, Pesticide Residues, Athens (2)

Recent advances in mass spectrometry have led to the development of new analytical methods applicable in food chemistry. Target and non-target mass spectrometry methods have been applied for the detection of known and unknown chemical contaminants, including pesticides and antibiotics residues in food. Mass spectrometry-based omics sciences have introduced high-throughput methods permitting the assessment of food authenticity and detection of adulteration. In this work, we focus mainly on mass spectrometry-based approaches such as metabolomics and lipidomics for studying small molecules in different food commodities (wine, honey, and olive oil) as well as LC-MS/MS and GC-MS/MS techniques for the detection of chemical contaminants at low concentration level.

ThP-94 / DEVELOPMENT OF INFANT FORMULA CERTIFIED REFERENCE MATERIAL FOR NUTRIENT ELEMENT ANALYSIS

Kyoung-Seok Lee - Jung Ki Suh - Sung Woo Heo - Sook Heun Kim - Youngran Lim - Myung Sub Han - Myung Chul Lim - Hyung Sik Min - Euijin Hwang - Yong-Hyeon Yim

KRISS, Center for Analytical Chemistry, Daejeon

Keywords: Certified reference material, Infant formula, Nutrient elements, Isotope dilution inductively-couple plasma mass spectrometry

Introduction

Nutritional status of infants and children has great impacts on subsequent health and disease over an entire lifetime. Supplying appropriate nutrition by well-balanced infant formula is especially important

considering that it is often the major, or even sole, source of nutrients at particular growth stages of an infant. It is a basic requirement for balanced nutrition and reliable nutritional studies to provide correct information on nutrients in foods and food supplements based on accurate measurements of their contents. Therefore, it is also important to develop a certified reference material (CRM) for the determination of elemental contents in infant formula or similar matrix. For this purpose, KRIS CRM 108-02-004 has been developed. The mass fractions of eight nutrient elements were certified by isotope dilution inductively coupled plasma mass spectrometry (ID ICPMS) and those of four monoisotopic elements were assigned as information values from the consensus values of expert laboratories. The present CRM is expected to provide a robust means for quality control and analytical method validation in the determination of various elements in infant formula or similar products.

Methods

Raw material of infant formula obtained from a manufacturer was spiked with six hazardous elements, Al, As, Cd, Hg, Sn and Pb, considering future expansion of the present CRM for hazardous element analysis as well as nutrient elements developed in this study. The spiked material was made into a paste form and blended thoroughly to make homogenous mixture. Then, the paste was freeze-dried, pulverized into powders, sieved to obtain particles between 50 to 250 μm , homogenized using V-mixer, bottled in 15 g each, and finally sterilized using 15 kGy ^{60}Co γ -ray irradiation. Twelve bottles from the batch of CRM were systematically selected and used for certification of nutrient elements. For certification of Ca, Cl, Cu, Fe, Mg, Se and Zn, a primary ratio method, ID ICPMS, was used. Monoisotopic element, I, Mn, Na and P, were analyzed using conventional field laboratory methods. Two expert laboratories in food analysis also participated in the analysis and their consensus values were given as information values instead of certified values.

Results

Mass fractions of eight elements certified by ID ICPMS showed excellent between-bottle homogeneity with relative standard deviation between 0.2 and 1.1 %. Expanded uncertainties of certified values were 1-4 %. In the case of monoisotopic element, Mn, Na and P, consensus values of KRIS and the two expert laboratories were given as information values. Iodine mass fraction was determined by gravimetric standard addition ICPMS with internal standard.

Conclusions

An infant formula CRM for element analysis was developed. Eight nutrient elements were certified, while 4 monoisotopic elements were given as information values. With additional certification including six hazardous elements, it will be the infant formula CRM with comprehensive coverage of analyte elements.

Novel Aspect

A new CRM for nutrient element analysis in infant formula certified with 8 elements with information values for 4 monoisotopic elements.

ThP-95 / ANALYSIS OF PGG(PENTA-O-GALLOYL β -D-GLUCOSE) IN HEALTH FUNCTIONAL FOODS BY LC-MS/MS

Shin Yong-Woon - Oh Bo Young - Jeong Hee Sun - Oh Jae Myoung - Hu Soo Jung - Hong Jin Hwan - Oh Keum Soon

Nutrition and Functional Food Research Team, National Institute of Food and Drug Safety Evaluation

PGG (Penta-O-galloyl beta-D-glucose) is currently used as functional ingredient in health functional food. This ingredient was approved by Korean Ministry of Food and Drug Safety (MFDS) as functional ingredient with the health benefit of the blood sugar level rises after eating. In this study, an analytical method based

on High performance liquid chromatography (HPLC) with positive ion electrospray ionization (ESI) coupled to tandem mass spectrometry detection (LC-MS/MS) was developed for the determination of PGG in health functional foods. Performance characteristics such as selectivity, linearity, recovery, repeatability, and intermediate precision were examined. This validated method was successfully applied to the determination of PGG in commercial functional foods.

ThP-96 / LIPID FINGERPRINT TO DETERMINE FUSARIUM CONTAMINATION IN MELON

Davila Zampieri⁽¹⁾ - ***Rodolfo Dantas***⁽¹⁾ - ***Katia Roberta Belaz***⁽²⁾ - ***Marcos N. Eberlin***⁽³⁾ - ***Andréia Hansen Oster***⁽⁴⁾ - ***Ebenezer de Oliveira Silva***⁽⁴⁾ - ***Guilherme Julião Zocolo***⁽⁴⁾

Federal University of Ceara, Organic and Inorganic Department, Fortaleza (1) - Federal University of Uberlândia, Institute of Chemistry, Uberlândia (2) - ThoMSon Mass Spectrometry Laboratory, UNICAMP, Campinas (3) - EMBRAPA, Embrapa Agroindústria Tropical, Fortaleza (4)

Keywords: Lipid fingerprinting, biomarkers lipid, Fusarium equiseti, Cucumismelo L., ESI-MS.

Introduction

The cultivation of melon (*Cucumis melo* L.) has become a socio-economic activity of the Brazilian semi-arid. This crop suffers serious damage due to phytopathogens, such as fungi of the genus *Fusarium* which cause melon rot, a disease that leads to total crop loss [1]. The MS has a technique that provides fast, efficient and reliable results to lipid fingerprinting has been for applications like clinical pathogen rapid and identification taxonomic classification of microbial collections [2,3].

Methods

The *Fusarium equiseti* was cultured in PDA for 14 days at 30°C (8 replicates). The pellet and melons were cut in 1x1 cm pieces. Lipids were extracted using CHCl₃/MeOH, according to Bligh Dyer method [4]. ESI-(+)-MS spectra were obtained using Xevo G2-XS Qtof Mass Spectrometer, covering the mass range from *m/z* 500 to 1000 Da. Peak picking and alignment were performed using MassLynx software and Multivariate analysis by MarkerLink (Waters).

Results

Fusarium equiseti was isolated from contaminated melons in the state of Ceará, Brazil. The lipid analyzes were obtained by ESI-(+)-MS, the exact masses and a fragmentation pattern allowed the use of the LIPID MAPS database and the injections in replicates allowed to perform multivariate analyzes showed characteristic ions of *Fusarium equiseti*, which are not present in the uninfected melon and which are found in the melon when it is infected by the phytopathogen. Based on the analyzes of PCA and O-PLS, the ions *m/z* 782.5 and *m/z* 784.5, these ions were identified as phosphatidylcholines (PC). Also observed the ions *m/z* 736.6 and *m/z* 783.5, identified as glycerolipids (GL), in the LIPID MAPS database. In all the replicates of the fungus samples, these ions were characteristic and were not found in the melon samples.

Conclusions

ESI-MS fingerprinting of lipid was able to classify if the melon is contaminated or uncontaminated by *Fusarium equiseti*. The PCA and O-PSL analyzes for these lipids, the most evident discrimination was for *Fusarium equiseti*, where highlighted ions which were attributed as GL and PC. Our preliminary results showed the potential the MS as a tool able for screening of major biomolecules participating in host-pathogen relationships.

Novel Aspect

To use the high efficient screening as quality control for melons, detecting contamination in the initial stage and avoiding losses of the lots during transportation.

References

1. Murad, N. B. A. Nor, N. M. I. M., Shohaimi, S., Zainudin, N. A. I. M. *Journal Applied Microbiology*, 123, 1533-1546, 2017.
2. Zhang, J. I., Talaty, N., Costa, A. B., Xia, Y., Tao, A., Bell, R., Callahan, J. H., Cooks, R. G. *International Journal of Mass Spectrometry*, 301, 37-44, 2011.
3. Willers, C., Jensen van Rensburg, P. J., Claassens, S. *Journal Applied Microbiology*, 118, 1251-1263, 2015.
4. Bligh, E. G., Dyer, W. J. *Canadian Journal of Biochemistry and Physiology*, 37, 911-917, 1959.

ThP-97 / INVESTIGATION OF THE GLIADIN FRACTION COMPOSITION IN OLD AND MODERN WHEAT GENOTYPES

Vincenzo Cunsolo⁽¹⁾ - *Michele A. De Santis*⁽²⁾ - *Antonella Di Francesco*⁽¹⁾ - *Rosaria Saletti*⁽¹⁾ - *Marcella M. Giuliani*⁽²⁾ - *Zina Flagella*⁽²⁾ - *Salvatore Foti*⁽¹⁾

Department of Chemical Sciences, University of Catania, Catania, Italy (1) - Department of the Science of Agriculture, Food and Environment (SAFE), University of Foggia, Foggia, Italy (2)

Keywords: proteomics; old wheat genotypes; gluten; allergens; celiac disease

Introduction

Over the last years, old wheat varieties have gained increasing attention since they have been suggested to present health benefits when compared with modern cultivars of bread and durum wheat and are generally considered better tolerated, even if without any scientific evidence based data [1]. Modern wheats have beneficial properties in terms of yield compared with older varieties, but information about changes in gluten proteins are still scanty.

Methods

In order to investigate the effects of the breeding that took place during the 20th century on gluten composition, a modern and an old Italian durum wheat genotypes were studied by coupling 2D-PAGE of the gliadin fraction and a bottom-up approach consisting of a specifically designed enzymatic digestion (trypsin plus chymotrypsin), high resolution mass spectrometry analysis and database search using a manually curated database (GluPro V1.0) of gluten proteins [3].

Results

A previous comparison of a set of old and modern Italian durum wheat genotypes [2] showed that 2D-PAGE profiles of the gliadin fractions did not present significant differences in relation to the expression of α - and γ -type gliadins. Conversely, they evidenced a drastic decrease in the expression of ω -type gliadins in the modern genotypes, with the highest content of ω -gliadins observed in the old genotype Dauno III and the lowest in modern genotype Saragolla. Particularly, the two genotypes presented differences in both the number and position of ω -gliadin spots. Here we report the preliminary results about the MS analysis of these spots. Taking into account that the main difficulties in ω -gliadin structural analysis are related to the presence of repeated peptide blocks rich of glutamine and proline, and to the very low level of lysine and arginine residues, a mixture of trypsin and chymotrypsin was employed with the aim to obtain an efficient ω -gliadin spots digestion.

Conclusions

The dedicated enzymatic digestion provided a lot of medium-sized peptides, which were analysed by MS/MS. Because proteomic methods require comprehensive and well annotated sequence databases, which are lacking for gluten, more reliable identifications were obtained searching MS data in a manually curated database of gluten proteins. This approach allowed the identification of some ω -gliadins responsible for the differences observed in 2D-PAGE profiles.

Novel Aspect

Provision of a custom approach together with the use of a curated and appropriate database will improve the reliability of gluten protein identification by proteomic analysis.

References

1. Bordoni A, Danesi F, Di Nunzio M, Taccari A, Valli V. *Int J Food Sci Nutr* 68, 278-286 (2016).
2. De Santis MA, Giuliani MM, Giuzio L, De Vita P, Lovegrove L, Shewry PR, Flagella Z. *Europ J Agron* 87, 19-29 (2017).
3. Bromilow S, Gethings LA, Buckley M, Bromley M, Shewry PR, Langridge JI, Mills ENC. *J Proteomics* 163, 67–75 (2017).

ThP-98 / NON-TARGETED CHARACTERIZATION OF THE DEGRADATION PRODUCTS OF SOYBEAN OIL AFTER FRYING PROCESS BY HS-SPME-GC-MS AND LLE-LC-HRMS

Maw-Rong Lee - Chung-Yu Chen - Ya-Ting Tsai

National Chung Hsing University, Chemistry, Taichung

Keywords: non-targeted analysis, degradation products, soybean oil, thermal processing

Introduction

Soybean oil has a maximum total production currently in edible oils and mostly used in the world. Furthermore, a high temperature thermal processing, such as frying and stir-frying have been proceeded in traditional Chinese cooking methods. The degradation compounds of soybean oil after thermal processing might cause the damage of human health like cancers [1]. Hence, to develop an efficiency analytical method for identification of degradation compounds of soybean oil after high temperature processing is important in food safety and toxicology. In this study, HS-SPME-GC-MS and LLE-LC-HRMS methods were developed for characterization of components produced in soybean oil before and after fried. Multivariate analysis has been widely used in non-targeted analysis in recent years. Principal components analysis (PCA) was utilized to assist the identification of the biomarkers used for distinguishing the fried and un-fried soybean oil.

Methods

The manual SPME holder and three kinds of SPME fibers including 65 μ m Polydimethylsiloxane/Divinylbenzene (PDMS/DVB), 75 μ m Carboxen/Polydimethylsiloxane (CAR/PDMS) and 50/30 μ m Divinylbenzene/Carboxen/Polydimethylsiloxane (DVB/CAR/PDMS) were obtained from Supelco (Bellefonte, PA, USA). These coated fibers were conditioned according to manufacturer's recommendation before the first use. 4 mL of soybean oil sample was placed in a 15 mL brown sample vial then sealed with a PTFE septum cap. Subsequently, the oil sample was heated to 175°C in oil bath with stirring bar at 300 rpm for imitating the heating processing while unheated oil at 50°C. Optimum incubation time was 3 min, and then extracted via headspace for 30 min at 175°C. After 30 min, the SPME fiber was desorbed in the GC injector. The separation of volatile analytes were performed with a TSQ Quantum GC system (Thermo Scientific, San Jose, CA, USA). Compounds were separated on a DB-5MS capillary column (30 m \times 0.25 mm i.d., 0.25 μ m film thickness). The injector was set at 270°C for 5 min. Ultra-high purity helium was used as the carrier gas at flow rate of 1 mL/min. The column temperature program was set as following: initial temperature 50°C for 2 min, subsequently raised to 300°C at the rate of 15 °C/min, then held at 300°C for 5 min at splitless mode. The GC/MS transfer line temperature was set at 275°C. Mass spectrums were obtained with full scan and mass range was from m/z 50 to 550. Library search was operated on NIST MS search 2.0.

Results

In order to extract all volatile compounds in soybean oil, the oil samples were incubated at 175°C for 3 min.

After incubation, the volatile components were extracted by CAR/PDMS, DVB/CAR/PDMS and PDMS/DVB SPME fibers at 175°C for 30 min as heated oil. Un-heated oil was processed at 50°C. The adsorbed extracts were desorbed in GC linear under 270°C for 5 min then analyzed by GC-MS/MS. For comparing the volatile compounds in standard soybean oil, brand A and B soybean oil, statistic technology was used for handling the GC-MS/MS data. PCA is most widely multivariate statistical method used to highlight similarities and differences. In this study, the obviously difference between heated and un-heated soybean oil was found by using PCA under 95% confidence interval. The volcano plot was used to discover the potential biomarkers for distinguishing the heated and un-heated soybean oil. The potential biomarkers was characterized by NIST database matching in the end. The matching results showed 1-Dodecene, 1-Undecanol, 1-Tridecanol, Hexanal, Tetrahydrofuran and 1-Tetradecanol were potential biomarkers in standard soybean oil; Nonanal, Dodecane, Benzaldehyde, 1-Octene, 1-Tridecanol and 1-Undecanol were potential biomarkers in brand A; Trimethylene oxide, Trichloromethane, Acetic acid, Tetrahydrofuran and 1-Tridecanol were potential biomarkers in brand B. In these compounds, Tetrahydrofuran, Trichloromethane were possible carcinogens. The biomarkers found in this study could be used for differentiation between the fresh soybean oil and recovered soybean oil.

Conclusions

This study concentrates on finding the possible biomarkers between fried and un-fried soybean oil by using HS-SPME-GC/MS and LLE-LC-HRMS. Use three kinds of fibers to make the data completed. We used heated oil to imitate fried oil. The experimental results indicated that there are significant differences between heated and unheated soybean oil in standard and two brands of oil. PCA will be employed to assist the classification and identification of biomarkers between fried and un-fried soybean oil. According to PCA and volcano plots, we can find a possible biomarkers by NIST. In the end, we find some possible biomarkers and google for their carcinogenicity. The experimental progress and results could assist us to distinguish fresh oil or not.

Novel Aspect

HS-SPME-GC-MS/MS and LLE-LC-HRMS combined with PCA was developed and used to characterize the biomarkers for soybean oil adulterating.

References

1. Inchingolo R., Cardenia V., Rodriuez-Estrada M. T., Lipid Technol., 25, 59-61 (2013)

ThP-99 / METABOLITES FROM CAFESTOL - A COFFEE DITERPENE- BY ZEBRAFISH TANK MODEL AND LC-MS/MS

Claudia M. Rezende⁽¹⁾ - *Cyrus V. Andriolo*⁽¹⁾ - *Henrique M. G. Pereira*⁽²⁾ - *Vinicius F. Sardela*⁽²⁾ - *Fabio J. M. Novaes*⁽¹⁾

Federal University of Rio De Janeiro, Chemistry Institute, Rio De Janeiro, Brazil (1) - Federal University of Rio De Janeiro, LBCD-Ladetec-Chemistry Institute, Brazil (2)

Keywords: mass spectrometry, metabolism, coffee diterpene, cafestol, in silico

Introduction

Coffee bean has a kaurane diterpene, cafestol [1], that acts as anti-inflammatory and anticancer (among others) [2], although being hypercholesterolemic [3]. It remains partially intact in the intestinal tract, and little is known about its metabolites. To understand what we really drink in the coffee beverage and its impact on health, this work was developed in an in vivo Zebrafish (*Danio rerio*) water tank model [4] with LC-MS/MS looking for cafestol metabolites.

Methods

The study was performed in 3 aquariums, 18 zebrafishes/aquarium and cafestol added at 1 mg ml⁻¹. The excretion water was submitted to SPE and then investigated in a Thermo Scientific Q Exactive Plus Orbitrap LC-MS/MS. The metabolites were characterized by at least 3 precursor ions. To corroborate the results and help rationalizing the fragmentation pathways observed, we compared the metabolites found with in silico analyses.

Results

Cafestol has a furan moiety linked to the ring A, a -CH₂OH (C17) and a tertiary OH both linked to C-16 at ring D of the diterpene moiety. All typical metabolic transformations were searched in the excretion water which, by tandem analyses, showed five cafestol metabolites with logical fragmentations. The most intense was proposed to be a partial epoxidation of the furan with a hydroxy group at C2 (*m/z* 349.20095). In other ion metabolites, we propose the intact furan with an aldehyde on C17 (*m/z* 337.17797); a -COOH at C17 (*m/z* 331.19038); a -OH at C2 besides -COOH at C17 (*m/z* 347.18585); and acetylation at -CH₂OH (*m/z* 359.22224). These showed that the preferential sites of metabolization were at carbons 2, 17, 18 and 19, compatible with in silico analyses with SMARTCyp and XenoSite. Although epoxidation was theoretically predicted, hydroxylation at C2 was not in the same metabolite.

Conclusions

Four metabolites observed in the water tank model of Zebrafish associated to LC-MS/MS are correlated to cytochrome P450 and one is a phase II metabolite. All transformations occurred at C- 2, C-17, C-18 or C-19 of the diterpene moiety, in agreement with in silico analyses.

Novel Aspect

This work contributes to a better understand of hot coffee beverage metabolites and helped on the fragmentation pathways of a complex structure like the diterpene cafestol, as also its metabolites.

References

1. Novaes, F. J. M. 2017. (DSc thesis). Rio de Janeiro. Federal University of Rio de Janeiro, Brazil.
2. George, S. E. et al. Critical Reviews in Food Science and Nutrition, 48, 464-486 (2008).
3. Urgert, R. et al. Journal of Agriculture and Food Chemistry, 43, 2167–2172 (1995).
4. Howe, K. et al. Nature, 496, 498 (2013).

ThPS-S06 - Lipidomics

Chairs: Gavin Reid, Donatella Caruso

ThP-100 / NOVEL SKIN LIPID BIOMARKERS ASSOCIATED WITH ATOPIC DERMATITIS

Evgeny Berdyshev⁽¹⁾ - ***Elena Goleva***⁽²⁾ - ***Irina Bronova***⁽¹⁾ - ***John Jung***⁽²⁾ - ***Max A Seibold***⁽²⁾ - ***Donald Y. M. Leung***⁽²⁾

National Jewish Health, Department of Medicine, Denver (1) - National Jewish Health, Department of Pediatrics, Denver (2)

Keywords: Atopic dermatitis; skin tape strips; sphingolipids; lysophosphatidylcholines

Introduction

Allergic skin diseases severely impair the quality of life of affected patients. Atopic Dermatitis (AD) is linked to a damaged skin barrier function that can precede clinical AD manifestations by many months. There are

no effective tools to predict the development of AD or to characterize the efficacy of applied therapy. Therefore, novel approaches to advance mechanistic understanding of AD as well as clinical management of AD and related skin diseases are needed.

Methods

A novel protocol for simultaneous assessment of lipids and filaggrin breakdown products in skin stratum corneum using tape stripping and liquid chromatography tandem mass spectrometry was employed. Targeted analyses with scheduled MRM algorithm were employed for precise quantitation of analytes. The hallmark of our approach is in the ability to normalize targeted mass spectrometric data per sample total protein content that was never done before.

Results

Lesional and non-lesional stratum corneum of 30 AD subjects was compared with 25 healthy individuals. In addition to the known redistribution of ceramide molecular species in AD, our study provided novel findings. Thus, the levels of sphingomyelins, and especially those with short chain fatty acids, increased several fold in atopic lesional skin. Most strikingly, lysophosphatidylcholines demonstrated a profound redistribution towards the prevalence of short chain species not only in lesional but also in non-lesional skin of AD subjects. RNA seq analysis performed on stratum corneum from same subjects revealed the downregulation of the expression of fatty acid elongates ELOVL3 and ELOVL6 that may explain the observed global redistribution of the long chain and short chain lipids. The modeling of the effect of type 2 cytokines on sphingolipids using human differentiated keratinocytes in vitro along with stable isotope pulse labeling confirmed the ability of type 2 cytokines IL-4 and IL-13 to provoke similar changes in sphingolipids.

Conclusions

Our ongoing studies further validate the power of the developed methodology and the value of lysophosphatidylcholines as biomarkers of atopic disease. Furthermore, skin tape stripping in conjunction with our unique processing protocol and targeted mass spectrometric methodology opens novel possibilities in clinical research by unveiling novel biomarkers of skin diseases and by allowing clinical studies in the skin using minimally invasive methodology.

Novel Aspect

Previously unrecognized association of sphingomyelins and lysophosphatidylcholines in skin stratum corneum with atopic dermatitis.

ThP-101 / QUANTUM CHEMICAL MASS SPECTROMETRY AS A TOOL FOR THE PREDICTION OF MASS SPECTRA AND THE IDENTIFICATION OF FRAGMENTATION PATHWAYS - PEPTIDES AND LIPIDS

Frank Blockhuys - Julie Cautereels

University of Antwerp, Department of Chemistry, Antwerp

Keywords: Quantum Chemical Mass Spectrometry, Tandem Mass Spectrometry, Density Functional Theory, Tripeptides, Lipids

Introduction

A new computational tool for the prediction of mass spectra based on quantum chemical calculations is being developed, called Quantum Chemical Mass Spectrometry for Materials Science (QCMS2). It was benchmarked by predicting the EI fragmentation pathways of a number of organics with a variety of functionalities. The main features in the mass spectra were correctly reproduced and new fragmentation routes proposed by QCMS2 were confirmed by MS/MS measurements [1].

Methods

QCMS2 is based on DFT/B3LYP/6-311+G* calculations and selects fragmentation pathways based on the relative values of bond orders, and reaction (for bond cleavages) and transition-state energies (for rearrangements). Known peptide-specific fragmentation mechanisms such as the Mobile Proton Model for His, confirmed by quantum chemical calculations, have been implemented into the method [2].

Results

QCMS2 has been applied to predict the fragmentations of tripeptides in ESI/MALDI CID MS, focusing on the influence of inter-side-chain (ISC) interactions on initial protonation and subsequent fragmentation. The fragmentations of non-cyclic tripeptides consisting of His as the central amino acid (AA) have been studied and a number of observations illustrate the importance of strong ISC interactions between the AAs; as a result, QCMS2 outperforms common tools such as PEAKS and Prospector in the assignment of tripeptide mass spectra [3].

Given the more recent interest in lipidomics QCMS2 is currently being used to predict the fragmentations of lipids. Common tools such as LipidXplorer and LipidBlast generally perform poorly due to the limitations of the models (e.g., the absence of McLafferty rearrangements). The results obtained for compounds such as sphingosylphosphorylcholine will be presented: fragmentation pathways in both positive- and negative-ion mode ESI/MALDI CID MS have been calculated and compared to the experimental spectra.

Conclusions

QCMS2 reproduces the main features in the mass spectra of tripeptides, *i.e.*, the traditional backbone cleavages resulting in a-, b₂- and y₁-ions, but, significantly, also predicts fragments resulting from ISC interactions. Likewise, it provides superior assignments of lipid spectra by taking a considerably larger number of fragmentation pathways into account. QCMS2 surpasses common tools in terms of insight into fragmentation and assignment of spectra.

Novel Aspect

QCMS2 assigns up to twice the number of signals in peptide/lipid spectra assigned by common mass spectral assignment software, and generates detailed fragment structures.

References

1. Cautereels J., Claeys M., Geldof D., Blockhuys F., J. Mass Spectrom., 51, 256-287 (2016).
2. Cautereels J., Blockhuys F., J. Am. Soc. Mass Spectrom., 28, 1227-1235 (2017).
3. Cautereels J., Van Hee N., Chatterjee S., Lemièrre F., Blockhuys F., Chem. Eur. J., submitted.

ThP-102 / A NEW UHPLC/MS-QTOF METHOD TO EVALUATE EFFECTS OF ENVIRONMENTAL METABOLIC DISRUPTORS ON HEPATIC LIPIDOMIC PROFILE

Fabrizia Carli⁽¹⁾ - **Melania Gaggini**⁽¹⁾ - **Veronica Della Latta**⁽¹⁾ - **Sabrina Tait**⁽²⁾ - **Francesca Maranghi**⁽²⁾ - **Roberta Tassinari**⁽²⁾ - **Cinzia La Rocca**⁽²⁾ - **Amalia Gastaldelli**⁽¹⁾

National Research Council (CNR), Institute of Clinical Physiology, Pisa (1) - Istituto Superiore di Sanità, Toxicology Unit, Roma (2)

Keywords: Ultra-High Pressure Liquid Chromatography Quadrupole Time of Flight (UHPLC/MS-QTOF), di-2-ethylhexyl phthalate (DEHP), hepatic lipidomic profile

Introduction

DEHP is a phthalate used in many consumer products. It is metabolized in the liver to MEHP, MEOHP and MEHHP. These compounds are known to alter organ lipid metabolism acting as metabolic disruptors[1]. Since the liver is a key organ in whole body metabolism and waste handling, the aim was to quantify DEHP metabolites in liver by UHPLC/MS-QTOF and relate them to alteration in organ lipid profile[2,3].

Methods

We evaluated the effect of 28 days exposure to DEHP on hepatic lipidomic profile. We studied 3 groups of rats (n=4 each, 2F/2M) exposed to low, medium and high dose vs non exposed. Lipids and DEHP metabolites were extracted from liver tissues (10-20 mg) and analyzed with ZORBAX Eclipse Plus C18 and ZORBAX SB Phenyl 2.1x100 mm 1.8 µm columns respectively. Concentrations were quantified using labeled internal standards added before the preparation.

Results

Hepatic concentrations of unconjugated DEHP metabolites (that reflects exposition) were increased after exposure. Male rats had significantly higher hepatic accumulation of both DEHP metabolites and lipids. Different lipid composition was found in rats exposed to DEHP vs non exposed: Ceramide CER(18:1/24:0) and CER(18:0:24:1) known to be involved in apoptotic pathway, were significantly higher ($p < 0.05$) in exposed rats indicating tendency to hepatic damage, while CER(18:2/20:0) was lower than in non-exposed rats ($p < 0.05$). These ceramides were synthesized more in proportion to total hepatic accumulation of DEHP metabolites. After DEHP exposure we also observed decreased hepatic phosphocholine PC(32:2) and phosphoethanolamine PE(38:2) ($p < 0.05$) and increased saturated triacylglycerols (TAGs) and diacylglycerols (DAGs), compounds known to be lipotoxic.

Conclusions

Rats exposed to DEHP showed a different hepatic lipid composition compared control rats, in particular in lipids involved in lipotoxicity and apoptosis.

Novel Aspect

We set up a new method by UHPLC/MS-QTOF to quantify DEHP metabolites and evaluate their metabolic effects from small amounts of livers of rats exposed to different doses of DEHP.

References

1. Frederiksen H., Skakkebaek Niels E., Andersson AM., Mol. Nutr. Food Res., 51, 899-911 (2007).
2. Rusyn I., Peters JM., Cunningham ML., Crit Rev Toxicol., 36, 459-479 (2006).
3. Saponaro C., Gaggini M., Carli F., Gastaldelli A., Nutrients, 7, 9453-9474 (2015)

ThP-103 / A QUANTUM CHEMICAL TOOL FOR THE IDENTIFICATION OF FRAGMENTATION PATHWAYS: ADDITIONAL INSIGHT INTO LIPID-SPECIFIC FRAGMENTATION MECHANISMS

Julie Cautereels - Frank Blockhuys

University of Antwerp, Department of Chemistry, Antwerp, Belgium

Keywords: Quantum Chemical Mass Spectrometry (QCMS2), Density Functional Theory (DFT), Tandem Mass Spectrometry, Tripeptides, Lipids

Introduction

A new computational tool for the prediction of mass spectra based on quantum chemical calculations is being developed, called Quantum Chemical Mass Spectrometry for Materials Science (QCMS2). It was benchmarked by predicting the EI fragmentation pathways of a number of organics with a variety of functionalities. The main features in the mass spectra were correctly reproduced and new fragmentation routes proposed by QCMS2 were confirmed by MS/MS experiments [1].

Methods

QCMS2 is based on DFT/B3LYP/6-311+G* calculations and selects fragmentation pathways based on the relative values of bond order, reaction energies (for bond cleavages) and transition state energies (for rearrangements). Know peptide-specific fragmentation mechanisms such as the Mobile Proton Model

(MPM) for His, confirmed by quantum chemical calculations, have been implemented into the method and the code omits the breaking of irrelevant bonds [2].

Results

QCMS2 was applied to predict the fragmentation patterns of tripeptides focusing on the fragmentation of the side chain and the influence of inter-side-chain (ISC) interactions on the fragmentation. Therefore, it is able to predict not only the traditional backbone cleavages but also side chain fragmentations and, as a result, convincingly outperforms common tools such as PEAKS [3].

Currently, QCMS2 is being used in the prediction of the fragmentations of lipids. Since the interest in lipidomics is more recent, there is not much known about their fragmentation and, therefore, the shortcomings of their methods are more pronounced than with peptides: none of the available methods can be used for all eight lipid classes and major fragmentation pathways such as McLafferty rearrangements are missing. The results of a number of sphingolipids including sphingosylphosphorylcholine will be presented: fragmentation pathways in both positive- and negative-ion mode ESI CID MS/MS have been calculated and compared to the experimental spectra.

Conclusions

QCMS2 correctly reproduces the main features in the mass spectra of lipids. Furthermore, QCMS2 associates, in contrast to literature, multiple(new) fragmentation pathways with the most prominent fragment ion, m/z 184 corresponding to phosphocholine in positive ion mode. In addition, QCMS2 predicts considerably more fragments in negative ion mode. Therefore, QCMS2 represents a considerable leap forward with respect to the available tools.

Novel Aspect

QCMS2 assigns a significantly higher percentage of signals in lipid mass spectra than the common tools, and provides detailed insight into the fragmentation behaviour.

References

1. Cautereels J., Claeys M., Geldof D., Blockhuys F., J. Mass Spectrom., 51, 602 (2016).
2. Cautereels J., Blockhuys F., J. Am. Soc. Mass Spectrom., 28, 1227 (2017).
3. Cautereels J., Van Hee N., Chatterjee S., Lemière F., Blockhuys F., submitted to Chem. Eur. J. (2018).

ThP-104 / ROBUST AND SENSITIVE LC-MS/MS BASED PLASMA LIPID PROFILING ON A THERMO SCIENTIFIC™ Q EXACTIVE™ HF-X MASS SPECTROMETER

Tabiwang N. Array ⁽¹⁾ – **Elena Sokol** ⁽²⁾ – **Angela Criscuolo** ⁽¹⁾ – **Clairey Daully** ⁽³⁾ – **Alexander Harder** ⁽¹⁾

Thermo Fisher Scientific, CMD, Bremen (1) - Thermo Fischer Scientific, CMD, Hemel Hempstead (2) - Thermo Fisher Scientific, CMD, Villebon sur Yvette (3)

Keywords: Lipidomics, Q Exactive™ HF-X, Profile Analysis

Introduction

Lipid profiling can be compromised by their extraction from biological samples, chromatographic separation, ionization and detection by mass spectrometry. The Thermo Scientific™ Q Exactive™ HF-X is a new Orbitrap based mass spectrometer with increased scan speed and sensitivity, which should both benefit lipid profiling experiments. In this study, we optimized a complete workflow for plasma lipidomics with improved sensitivity and quantitative performances.

Methods

Human plasma samples were spiked with known concentrations of SPLASH™ Standard (Avanti) and

extracted with varying sample amounts. Lipids were separated on a 15 cm Accucore™ C30 column and analyzed on the Q Exactive HF-X. Several instrument parameters were tested, including detection mass range, ion source conditions, ion transfer and data dependent parameters. LipidSearch™ software was used for identification and quantification of molecular lipid species.

Results

Different plasma extraction volumes were tested in order to maximize the number of lipid species which could be detected and quantified with good precision. This evaluation was conducted by extracting the ion signals for the 14 spiked isotopic labeled lipid standards representing all major lipid classes present in plasma. Replicate extractions confirmed the reproducibility of the method with coefficient of variations below 20% for all standards. Further optimizations included the limitation of in-source fragmentation of labile lipids and improvement of MS1 detection across the studied lipid classes. Calibration curves were generated to determine the linearity of detection and limits of quantification. Finally, data dependent parameters were evaluated in order to maximize the number of high quality MS2 spectra generated in view of the higher scan rate and sensitivity of the instrument.

Conclusions

This method optimization resulted in the confident identification and quantification of several hundreds of lipids from a human plasma sample. For Research Use Only.

Novel Aspect

Development of a robust and sensitive mass spectrometry based workflow for large-scale plasma lipidomics profiling studies.

ThP-105 / VERIFIED HILIC LC-MS/MS ASSAY FOR HIGH-THROUGHPUT TARGETED LIPDOMICS ANALYSIS

Cyrus Papan⁽¹⁾ - *Santosh Kapil Kumar Gorti*⁽²⁾ - *Lijuan Fu*⁽³⁾ - *Baljit K. Ubhi*⁽³⁾ - *Mackenzie J Pearson*⁽³⁾ - *Lei Xiong*⁽³⁾ - *Paul Baker*⁽³⁾ - *Jörg Dojahn*⁽¹⁾

SCIEX Germany GmbH, Sciex, Darmstadt (1) - SCIEX USA, Sciex, Framingham, MA (2) - SCIEX USA, Sciex, Redwood City, CA (3)

Introduction

There are about 150,000 different lipid molecular species present across the biological spectrum. 'Shotgun' lipidomics is an established non-targeted approach for broad-based lipidomic analysis but it can suffer from inherent ion suppression and ambiguous identification. Herein we introduce a targeted lipidomics workflow with the combination of HILIC LC separation and multiple reaction monitoring (MRM) mass acquisition. HILIC LC chemistry provides reproducible isomer separation based on lipid classes. The Scheduled MRM™ type data acquisition provides high sensitivity and wide coverage (over 1200 lipid species) for lipids screening and relative quantitation with 24 min run time. The method was extensively verified and is a robust solution for a targeted lipidomics workflow.

Methods

Standards and samples Lipids standards (SPLASH heavy isotopic labeled standards, 17:1 standards and regular light standards) and bovine heart extract were purchased from Avanti and SCIEX. They were diluted with 50%:50%/MeOH: dichloromethane containing 10mM ammonium acetate before injection. LC-MS Bovine heart extract samples were subjected to LC-MS/MS analysis with a QTRAP® 6500+ mass spectrometer coupled with Exion HPLC system. A Waters XBridge Amide 3.5 μm, 4,6 x 150 mm column was used for HILIC based LC separation. A Scheduled MRM™ method was implemented to provide screening and relative quantification of over 1200 lipids with 24 min run time.

Preliminary Data

Isomer interference minimization among different lipid classes is known as one of major challenges for LC-MS/MS method development for lipidomics analysis. Very commonly, one lipid molecule can have multiple isomers in other lipid classes sharing the identical precursor and product masses. This cross-class crosstalk could be troublesome because mass spectrometers cannot differentiate these isomers based on MRM transitions. To overcome this issue, a HILIC chromatographic method is implemented to provide good LC separation based on lipid classes. In order to confirm separation efficiency, lipids standards as natural lipids extracts distinguished by lipid class (one class per standard) were injected individually to confirm that there was no isomer crosstalk among different lipid classes. Retention time reproducibility improvement retention time on a HILIC chromatography is very sensitive to changes of LC condition. A minor pH difference or organic composition change in the mobile phase might induce retention time shifts of analytes. To minimize retention time drift, extensive evaluation on mobile phase preparation was performed. A standard LC preparation protocol was developed to achieve optimal retention time reproducibility, which is especially critical for the Scheduled MRM™ method. Internal Standard Selection Lipid standards from 19 different classes, which are either heavy isotopic labeled lipids or odd chain lipids, served as internal standard. This method provided extensive lipid class coverages including, CE, CER, DCER, HCEr, LCER, TAG, DAG, MAG, LPC, PC, LPE, PE, LPG, PG, LPI, PI, LPS and PS.

Novel Aspect

HILIC-MS/MS by QTRAP® system provides a reliable solution for targeted lipidomics workflow with maximum sensitivity, wide coverage and high throughput.

ThP-106 / STRUCTURE ELUCIDATION OF GRAM-NEGATIVE MICROBIAL GLYCOLIPIDS BY ENERGY-RESOLVED TANDEM MASS SPECTROMETRY

Agnes Dornyei⁽¹⁾ - Viktor Sandor⁽²⁾ - Aniko Kilar⁽²⁾ - Nikolett Eros⁽¹⁾ - Bela Kocsis⁽³⁾ - Ferenc Kilar⁽²⁾

University of Pecs, Department of Analytical and Environmental Chemistry and Szentagotthai Research Centre, Pecs (1) - University of Pecs, Institute of Bioanalysis, Pecs (2) - University of Pecs, Department of Medical Microbiology and Immunology, Pecs (3)

Keywords: bacterial lipid A, structure elucidation, fragmentation pathway, energy-resolved tandem mass spectrometry

Introduction

Gram-negative bacteria are responsible for a variety of foodborne diseases and healthcare-associated infections. The outer leaflet of their outer membrane is composed of lipopolysaccharides (LPSs) of diverse structure and composition. The beneficial (immune-stimulatory) and/or adverse (*e.g.*, pro-inflammatory) effects of these molecules released from the bacteria depend mainly on the structure of the lipid A, the lipophilic anchor constituent of the LPS.

Methods

The HPLC-MS/MS measurements were performed on an Infinity 1290 UHPLC system coupled to a 6530 Accurate Mass Q-TOF mass spectrometer (Agilent) equipped with a Jet Stream ESI ion source. The separation was achieved using a core-shell Kinetex™ C18 column. The collision energy dependence of the low-energy CID fragmentation pathways was followed both in negative- and positive-ion modes. Lipid A samples obtained from *Escherichia coli* strains were used as “training sets”.

Results

Recently, we have established a new RP-HPLC method combined with ESI-QTOF MS/MS for the simultaneous determination and structural characterization of complex and heterogeneous lipid A samples [1-3]. We have also investigated the effects of the application of different eluent additives on the

separation and on the ionization. The fragmentation pathways of the adduct ions and of the protonated lipid-A molecules in the positive-ion mode and those of the deprotonated molecules in the negative-ion mode were explored and compared. Phosphorylation site-specific diagnostic fragment ions were identified in the positive-ion mode. By a comprehensive analysis of the fragmentation pathways of numerous lipid A variants, a full structure elucidation method was proposed.

Conclusions

The overall structural analysis of lipid A samples is of high importance, because this knowledge is essential to describe any structure-function relationship related to Gram-negative bacterial endotoxins. Energy-resolved tandem mass spectrometry proved to be a valuable method to mapping trends in fragmentation relative to the collision energy applied during the activation.

Novel Aspect

We present an energy-resolved tandem mass spectrometry based strategy to characterize the microbial membrane lipid A structures.

References

1. Sandor V., Kilar A., Kilar F., Kocsis B., Dornyei A., Journal of Mass Spectrometry, 51, 615 (2016).
2. Sandor V., Dornyei A., Makszin L., Kilar F., Peterfi Z., Kocsis B., Kilar A., Journal of Mass Spectrometry, 51, 1043 (2016)
3. Sandor V., Kilar A., Kilar F., Kocsis B., Dornyei A., Journal of Mass Spectrometry, 53, 146 (2018).

Acknowledgements

The research was supported by the ÚNKP-17-4-III New National Excellence Program of the Ministry of Human Capacities, and the grant NKFIH K-125275.

ThP-107 / DEVELOPMENT OF A COMPREHENSIVE PANEL OF EICOSANOIDS AND RELATED FATTY ACID METABOLITES FOR HIGH-SENSITIVITY AND QUANTITATIVE LC-MS/MS ANALYSIS

Yuka Fujito⁽¹⁾ - *Masaki Yamada*⁽¹⁾ - *Huan Lin*⁽¹⁾ - *Setsuko Fujita*⁽²⁾

Shimadzu Corporation, Global Application Development Center, Kyoto (1) - Ono Pharmaceutical Co., LTD., Discovery Technology Research Laboratories, Osaka (2)

Keywords: Eicosanoid, lipid mediator, omega-3 fatty acid, MRM, ionization polarity switching

Introduction

Eicosanoids or omega-3 fatty acids metabolites are attractive targets to understand physiological roles. Comprehensive and quantitative monitoring methods using LC/MS has been developed in last decades. Here we report a novel method consisted of 300 MRM transitions packed in a 20 minutes' chromatographic condition for 200 targets. From 30 µL human serum, 75 targets were assigned by the consistence of retention time or intensity ratio with authentic standards.

Methods

Authentic standards were purchased from Cayman Chemical (Ann Arbor, MI). Human serum and blood were obtained from Sigma-Aldrich Co. LLC. and from Bioreclamation IVT, US, respectively. STRATA-X (Phenomenex, Torrance, CA) cartridge were used for solid phase extraction. An LC/MS system consisting of Nexera UHPLC system and LCMS-8060 (Shimadzu Corp.) was used. A Kinetex C8 column (2.1 x 150 mm, 2.6 µm) was used for target separation.

Results

A new method consisted of over 300 MRM transitions for 186 eicosanoids and related fatty acid metabolites and of 14 deuterium labeled analogs as internal standards. Expected retention times were confirmed by analyzing authentic standards with a 20 minutes' chromatographic condition using 0.1% acetic acid and acetonitrile as elution buffers. More than one MRM transitions were applied for over 100 targets. The most of MRM transitions were negative ion mode without 27 transition with positive ion mode. When 5 msec dwell time and 1 msec pause time were set, lower limit of quantifications were sub-pg for the majority of targets. The quantitative profiling of human serum resulted that 75 targets were assigned by the consistence of retention time within 3 seconds or that of the intensity ratio detected in multiple MRMs with corresponding authentic standards. Major eicosanoids PGE2 and PGD2 were not detected. An arachidonic acid metabolite 11-HETE was observed with the highest intensity close to saturation of the detector.

Conclusions

We have developed a new LC/MS method for comprehensive quantification of over 180 eicosanoids and related fatty acid metabolites. The method was applied to 30 µL of human serum sample to assign 75 targets including 30 arachidonic acid metabolites, 13 DHA metabolites, 9 EPA metabolites. The method will be detectable over 100 targets which covers wide range of fatty acid metabolites. We believe the method is useful to study physiological functions or biomarker discovery.

Novel Aspect

We successfully developed a comprehensive MRM panel of 186 eicosanoids and related fatty acid metabolites for high-sensitivity and quantitative LC-MS/MS analysis.

References

1. Kita, Y. ; Takahashi, T. ; Uozumi, N. ; Shimizu, T. Analytical biochemistry, 342, 134-143 (2005).
2. Yamada, M. ; Kita, Y. ; Kohira, T. ; Yoshida, K. ; Hamano, F. ; Tokuoka, M. S. ; Shimizu, T. J. Chromatography B, 995, 74-84 (2015).

ThP-108 / RAPID TARGETED LC-MS/MS ASSAY FOR ESKAPE PATHOGEN IDENTIFICATION DIRECT FROM BIOLOGICAL FLUIDS

David Goodlett⁽¹⁾ - ***Tao Liang***⁽¹⁾ - ***Benjamin Oyler***⁽²⁾ - ***Courtney Chandler***⁽³⁾ - ***Sung Hwan Yoon***⁽³⁾ - ***Robert Ernst***⁽³⁾

UMB, Pharmacy, Baltimore (1) - UMB, IMET, Baltimore (2) - UMB, Dentistry, Baltimore (3)

Keywords: pathogen identification, rapid, polymicrobial, antibiotic

Introduction

Microbial identification relies on bacterial culture followed by biochemical testing, which generally takes overnight. Protein analysis by MALDI-TOF MS has become a prominent player in clinical microbial diagnosis. Identification is nearly instantaneous, but requires overnight growth to produce the required pure colony. We recently published (1,2) a method for microbial identification direct from biofluid without the need for culture.

Methods

Bacterial membrane glycolipids were extracted as published (1) and reconstituted in chloroform/methanol/water (3: 1.5: 0.25, v/v/v) for LC-MS/MS analysis on a SYNAPT G2 (Waters). Signature ions from each ESKAPE species were selected and fragmented. Suitable transitions were determined for each ESKAPE pathogens. The existing MALDI library (1) was used to confirm LC-MS results after deconvolution and

manual inspection used for polymicrobial analysis.

Results

Previously we published use of microbial lipids for identification of bacteria and fungi direct from specimen using the Bruker Biotyper platform (1, 2). Here we present development of a pseudo-MRM assay for the ESKAPE pathogens (1) on a SYNAPT MS. A new extraction protocol that takes less than one hour is also introduced. Individual glycolipid preparations were analyzed by direct infusion to generate mass spectra for each ESKAPE species. Signature lipid ions, which are those unique to a given microbe, were selected for tandem MS fragmentation and their unique transitions, attributable to loss of various fatty acids, noted. For example, precursor ions at m/z 1744, 1824, and 977 ($z=2$) were found to be unique to *Enterobacter cloacae*. After successfully generating tandem mass spectra for each individual ESKAPE pathogen, mixed samples containing multiple pathogens were tested to demonstrate the method could distinguish and detect each organism in a simulated polymicrobial sample and where appropriate antibiotic resistance (3).

Conclusions

ESKAPE pathogens, *Enterococcus faecium*, *Staphylococcus aureus*, *Klebsiella pneumoniae*, *Acinetobacter baumannii*, *Pseudomonas aeruginosa*, and *Enterobacter* spp., were targeted due to their high incidence of nosocomial infections and multidrug resistance. Here we present initial work to develop a rapid (< 1 hour) pseudo-MRM assay targeted at detection of these microbes that are sensitive and resistant to antibiotic therapy.

Novel Aspect

Rapid identification of ESKAPE pathogens by LC-MS/MS points toward a targeted clinical assay that may be carried out on triple quadrupole systems rather than MALDI-TOF MS.

References

1. Leung LM, Fondrie WE, Doi Y, Johnson JK, Strickland DK, Ernst RK, Goodlett DR. *Scientific Reports* 25, 7, 6403 (2017).
2. Liu YY, Chandler CE, Leung LM, McElheny CL, Mettus RT, Shanks RMQ, Liu JH, Goodlett DR, Ernst RK, Doi Y. *Antimicrob. Agents Chemother.* 24, 61 (2017).
3. Leung LM, Cooper VS, Rasko DA, Guo Q, Pacey MP, McElheny CL, Mettus RT, Yoon SH, Goodlett DR, Ernst RK, Doi Y. *J. Antimicrob. Chemother.* 72, 3035-3042 (2017).

ThP-109 / BIOMARKER RESEARCH AND SHORT TIME CHROMATOGRAPHY: A UTOPIA?

Lisa Hahnfeld - Gerd Geisslinger - Nerea Ferreirós

Goethe University Frankfurt, Institute of Clinical Pharmacology, Frankfurt

Keywords: LC-HRMS, lipidomics, biomarker research

Introduction

In biomarker research short analysis times are preferred, which are limited by problems like matrix effects and coelution of isobaric compounds. For this reason, run-times often exceed 30 or even 50 minutes. Our aim is to implement an untargeted method for high-throughput screening, reducing drifts in instrumental response, retention time shifts and improving sample stability.

Methods

We developed a chromatographic method with a sub 2 micrometer C8 column and a binary gradient to cover the lipidome ranging from polar lipids to highly apolar triacylglycerides and cholesteryl esters in biological samples after liquid-liquid extraction. The chromatographic system was coupled with quadrupole-time-of-flight mass spectrometry for lipid identification.

Results

The obtained chromatographic resolution allowed the separation of isobaric lipids and reduced ion suppression effects. Combined with high resolution mass spectrometry we could identify lipids from over 15 different lipid classes. The resulting acquisition time was 17 minutes per sample and polarity, simplifying the measurement of large sample numbers and improving system stability. Hence we achieved a coefficient of variation (CV) below 15% for the peak area of 220 lipids in human plasma at repeated extractions, with 205 of these having a CV below 10%. Additionally, the average CV of the retention times was 0.1%. We successfully applied the method on plasma, serum and cell culture samples as well as several different mouse tissue types. For the latter, we could show the importance of the sampling location of mouse kidney, heart and liver for the reproducibility of the lipid profile.

Conclusions

The developed method can be applied for high-throughput lipidomics and biomarker discovery research in biological samples.

Novel Aspect

LC-HRMS analysis of biological samples in 17 minutes with low matrix effects, adequate peak resolution, lipid identification and method reproducibility.

ThP-110 / INVESTIGATION INTO THE BLEACHING OF VEGETABLE OILS USING HIGH RESOLUTION MASS SPECTROMETRY

Andrew Hambly - Jeroen Van Duijneveldt - Paul Gates

University of Bristol, School of Chemistry, Cantock's Cl, Bristol BS8 1TS, Bristol

Keywords: HPLC-MS, vegetable oils, bleaching, adsorbents

Introduction

The refining of lipid-based oils is important to the food and pharmaceutical industries in the production of food supplements and excipients. 'Bleaching' refers to one step in the process where the oils are mixed with adsorbents to remove contaminants such as; lipid peroxides, pigments and soaps. This project uses high resolution mass spectrometry to study the bleaching of vegetable oils.

Methods

Various techniques including reverse-phase UHPLC-MSn and MALDI MS were used to analyse crude and refined versions of sesame, olive and fish oils, as well as synthetic poly(ethylene glycol) (PEG) and polysorbate based products. A rapid method for processing spectra was built using software package Knime.

Results

Refined versions of vegetable oils saw lower levels of lipid oxidation products and pigments as predicted from drops in peroxide value and removal of colour after processing. Refined oils also showed lower levels of soaps, metals and phospholipids. Compounds seeing the greatest removal saw a drop in abundance by factor 20.

Surprisingly, refined oils also showed lower concentrations of diglycerides. This was suggested to be due to presence of the glyceryl hydroxyl group.

Conclusions

Combination of high resolution mass spectrometric methods with computer processing can provide a high throughput means of studying the removal of compounds during edible oil processing. Data confirmed removal of typical contaminants during oil processing. It also highlighted change of otherwise unknown

compounds during bleaching.

Novel Aspect

Use of high resolution mass spectrometric methods as a means of rapid study of bleaching processes in edible oil refinement.

ThP-111 / ANALYSIS OF A COMPLEX LIPID PROFILE IN HUMAN SEBUM BY USING GCXGC-HRTOFMS

Masahiro Hashimoto⁽¹⁾ - *Koji Okuda*⁽²⁾ - *A. John Dane*⁽²⁾ - *Robert B. Cody*⁽²⁾ - *Bram van der Meer*⁽³⁾

JEOL (EUROPE) SAS, European Application Group, Croissy-sur-Seine (1) - JEOL USA Inc., Application Group, Peabody, MA (2) - JEOL (Europe) B. V., Sales Executive Science and Measurement Instruments, Nieuw-Vennep (3)

Keywords: Lipid analysis, High resolution, High mass accuracy

Introduction

Sebum is a complex mixture of lipids that are typically measured by using GC-FID. However, these methods are often unable to separate the individual components from each other which makes it difficult to definitively identify each of them. Conversely, GCxGC-HRTOFMS is a powerful tool for identifying analytes in complex mixtures. In this work, we measured lipids in sebum by using GCxGC-HRTOFMS.

Methods

The sample was collected onto square aluminum sheets wiped across a human forehead. The lipid compounds were then extracted from the aluminum sheets using a 50/50 volume solution of methanol/dichloromethane with sonication. Next, the extraction solution was centrifuged and concentrated. Afterwards, the final solution was analyzed using GCxGC-HRTOFMS.

Results

A lot of lipid compounds such as free fatty acids, wax esters, diacylglycerols and triacylglycerols were observed, and these lipid compounds were separated by GCxGC perfectly. Obtained peaks were comprehensively identified using a combination of accurate mass measurements and library searches. As an example, the peaks detected at retention time 18.684 minute (1st column), 6.756 seconds (second column) was matched to oleic acid using a NIST library search. Furthermore, this mass spectrum showed a molecular ion that coincided with the calculated accurate mass of oleic acid (error of 1.4 mDa). The results for the other compounds present in the sebum sample will be presented in this work.

Conclusions

Low boiling point fatty acids up to high boiling triglycerides were observed in a single sebum measurement. The cholesterol esters and triacylglycerols, which are difficult to separate by 1D GC, showed high chromatographic separation by GCxGC-HRTOF. Additionally, the peaks were comprehensively identified by using a combination of accurate mass measurements and library searches.

Novel Aspect

Analysis of a complex lipid profile in human sebum by using GCxGC-HRTOFMS.

ThP-112 / APPLICABILITY OF MALDI-ORBITRAP-MS FOR THE EARLY DIAGNOSIS OF RENAL CELL CARCINOMA BASED ON THE LIPIDOMIC ANALYSIS OF BLOOD AND URINE SAMPLES

Robert Jirásko⁽¹⁾ - *Denise Wolrab*⁽¹⁾ - *Ivana Brabcová*⁽¹⁾ - *Michal Holčapek*⁽¹⁾ - *David Vrána*⁽²⁾ - *Vladimír Študent*⁽³⁾ - *Bohuslav Melichar*⁽²⁾

University of Pardubice, Department of Analytical Chemistry, Pardubice (1) - Palacký University, Department of Oncology, Olomouc (2) - Palacký University, Department of Urology, Olomouc (3)

Keywords: MALDI, Orbitrap, lipidomics, sulfatides, renal cell carcinoma

Introduction

Renal cell carcinoma (RCC) is the most common type of kidney cancer. The RCC incidence and the mortality rate have been increasing globally during the last decade. Up to date, there are still no recommended screening tests for early diagnosis. In our work, the matrix-assisted laser desorption/ionization coupled with Orbitrap mass spectrometry (MALDI-Orbitrap-MS) is used for the clinical study of patients with RCC and healthy volunteers.

Methods

All samples were spiked before the extraction with the appropriate internal standard mixture. Lipidomic extracts of human plasma were prepared by modified Folch method, while reversed-phase solid phase extraction was applied for urine samples. The LTQ Orbitrap instrument was operated in the negative-ion mode over a mass range m/z 400-2000 and the mass resolution was set to 100,000. 9-aminoacridine was used as MALDI matrix.

Results

First, lipid species are identified using high mass accuracy full scan and tandem mass spectra. Subsequently, the optimization and method validation for semiquantitative measurement is performed. The optimized methodology is finally used for the analysis of all studied body fluid samples. Each measurement is represented by one average MALDI-MS spectrum with thousands of m/z peaks. The automatic peak assignment is subsequently performed and particular m/z peaks of various intensities are matched with deprotonated molecules from a database created during the identification procedure using the Excel macro script for each sample. Multivariate data analysis methods, *i.e.*, unsupervised principal component analysis and supervised orthogonal partial least square discriminant analysis, are further used for the visualization of differences between studied sample groups.

Conclusions

Most upregulated and downregulated lipids in blood and urine samples of RCC patients are finally visualized using particular boxplots. Significant differences in lipid abundance between RCC patients and healthy volunteers proved potential of MALDI-Orbitrap-MS for early diagnostic of RCC, but further testing is still necessary. This work was supported by the Czech Science Foundation (project No. 18-12204S).

Novel Aspect

Applicability of MALDI-Orbitrap-MS for the early diagnosis of renal cell carcinoma based on the lipidomic analysis of blood and urine samples.

ThP-113 / CHANGES OF SPECIFIC SPHINGOMYELINS DEPENDING ON THE SEVERITY OF NONALCOHOLIC FATTY LIVER DISEASE

Youngae Jung ⁽¹⁾ - Jungju Seo ⁽²⁾ - Geum-Sook Hwang ⁽¹⁾

KBSI, Western Seoul Center, Seoul (1) - Scientific Instruments Reliability Assessment Center, Korea Basic Science Institute (2)

Keywords: Nonalcoholic fatty liver disease, lipidomics, sphingomyelin

Introduction

Nonalcoholic fatty liver disease (NAFLD) is one of common causes of chronic liver disease that involve hepatic accumulation of triglycerides. The clinical-histologic spectrum of NAFLD extends from a nonalcoholic fatty liver (NAFL) to nonalcoholic steatohepatitis (NASH). NASH is divided from NAFL according to the presence of ballooning and lobular inflammation [1]. Although NAFLD is closely correlated with obesity, 5-8% of lean subjects also have NAFLD [2].

Methods

Total 361 serum samples from patients with biopsy-proven NAFLD were classified in six groups on the basis of BMI (25 kg/m²) and liver histology (No NAFLD, NAFL and NASH). Lipidomics study for all serum samples were performed using ultra-performance liquid chromatography/quadrupole time-of-flight mass spectrometry (UPLC/QTOF MS).

Results

The diverse lipid metabolites identified in serum samples showed a characteristic alterations depending on the severity of NAFLD between non-obese and obese groups. Most of all, specific sphingomyelins (SMs) had distinct increase patterns according to the stage of steatosis, lobular inflammation, and ballooning except for fibrosis only in non-obese NAFLD. Moreover, the receiver operator characteristic (ROC) analysis for determining the progression degrees of NAFLD produced higher the area under the curve (AUC) values in non-obese group than in obese.

Conclusions

These results presented that disease mechanism is different between non-obese and obese NAFLD, and these specific SMs might be used as markers to predict the severity and/or progression of NAFLD.

Novel Aspect

Specific SMs were clearly associated with liver histology including steatosis, lobular inflammation, and ballooning in non-obese NAFLD patients.

References

1. Puri P., Baillie RA., Wiest MM., et al., *Hepatology*, 46, 1081-1090 (2007).
2. Feldman A., Eder SK., Felder TK., et al., *Am. J. Gastroenterol.*, 112, 102-110 (2017).

ThP-114 / A COMPREHENSIVE METHOD FOR TOTAL LIPID ANALYSIS IN BIOLOGICAL MATRICES BY HYDROPHILIC INTERACTION LIQUID CHROMATOGRAPHY ELECTROSPRAY IONIZATION MASS SPECTROMETRY HILIC/ESI-MS

Spiro Khoury - Stéphanie Cabaret - Elodie A. Y. Masson - Olivier Berdeaux

UMR CSGA - INRA - CNRS, University of Bourgogne Franche-Comté, Dijon

Keywords: Lipids, Extraction, liquid chromatography, Mass spectrometry, Ion suppression

Introduction

Lipid analyses are increasingly used for determining the development of several diseases. Because of their extremely heterogeneous structures, analysis of lipids is challenging. LC/ESI-MS, the major tool used to analyze complex lipidomes, is very sensitive towards the composition of lipid extract. Therefore, a complete analytical strategy is proposed in this work in order to efficiently extract, purify and analyze lipids in different biological matrices.

Methods

Lipids were extracted from rat brain and plasma samples with an optimized method using a mixture of

solvents of increasing polarity. The efficiency of this method was compared to other reference methods. Then, lipids of low relative abundance were purified by several methods (C18 solid-phase, phase partition [1] or Phree phospholipid removal columns). LC-MS (Triple quadrupole and Orbitrap FUSION) methods were developed for targeted and global analyses of lipids.

Results

The optimized monophasic extraction method enabled a fast and efficient extraction of all lipids from the matrices. A better recovery of lysophospholipid and ganglioside (GG) classes, whereas similar recovery of other phospholipids (PL), ceramides, fatty acids and sterol, was obtained compared to other extraction methods. The HILIC conditions allowed the separation of PL and GG classes, including acetylated forms, according to their increased polarity. Analyses conducted by LC-ESI-HRMS/MS enabled the identification of a large diversity of lipid molecular species. Quantification of lipids by MS tuned to specific acquisition modes showed high levels of phosphatidylcholine and phosphatidylethanolamine classes. Different purification methods of GG were evaluated in order to overcome some analytical challenges in LC/ESI-MS like ion suppression phenomenon [2] after co-elution with the more abundant PL of similar polarity. Therefore, an increase in quantification of all GG classes was observed after removing PL prior to analysis.

Conclusions

Polar and complex lipids in biological samples usually present very low extraction yield with the most popular Folch-type procedures. In this work, we propose a monophasic extraction method adapted to extract lipids of a large range of polarity. A method in tandem MS was optimized and coupled to an HILIC system to separate, characterize and quantify all lipids. In this respect, powerful purification tools were used to avoid ion suppression effect for ESI-MS analysis.

Novel Aspect

Our complete analytical strategy offers an efficient extraction of both polar and non-polar lipids from complex biological matrices followed by an accurate LC-MS analysis.

References

1. A. Y. Masson E et al., J Lipid Res, 56, 1821-35 (2015)
2. Khoury S et al., Anal Bioanal Chem, 408, 1453-65 (2016)

ThP-115 / GRAPHICAL IDENTIFICATION OF LIPIDS USING THREE-DIMENSIONAL KENDRICK MASS PLOTS

Ansgar Korf - Christian Vosse - Heiko Hayen

University of Münster, Institute of Inorganic and Analytical Chemistry, Münster

Keywords: Lipidomics, Kendrick mass defect, High-resolution mass spectrometry, Liquid chromatography

Introduction

Recent technical advances regarding liquid chromatography (LC) and HRMS enable the mapping of the lipidome of an organism with short data acquisition times. However, interpretation and evaluation of resulting multidimensional datasets are challenging and this is still the bottleneck regarding overall analysis times. Therefore, we present a novel adaption of Kendrick mass (KM) plot analysis.

Methods

Separation of lipids by their respective head group was achieved via hydrophilic interaction liquid chromatography (HILIC) coupled to a HRMS instrument (Q Exactive). The resulting LC/HRMS datasets were processed to a list of chromatographically separated features by applying an optimized MZmine 2 workflow. [1] All features were plotted in a 3DKM plot for graphical lipid identification.

Results

To demonstrate the novel adaption of KM plot analysis, occurring phospho- and glycolipid classes of the green alga *Chlamydomonas reinhardtii* have been separated using HILIC. Phospho- and glycolipid species have been identified by exact mass database matching and graphical 3DKM plot analysis according to [2]. LC-HRMS/MS datasets have been processed utilizing an optimized MZmine 2 workflow for feature list generation. Database matching was achieved by extending the already existing MZmine 2 module 'glycerophospholipid search' by adding missing phospho- and glycolipid classes. Furthermore, graphical analysis has been carried out by implementing a new visualization module in MZmine 2, which adds chromatographic characteristics, such as retention time, to a classic KM plot in form of a heat map. The developed visualization module is part of the MZmine 2 toolbox since version 2.31.

Conclusions

The use of the novel adaption of the KM plot has accelerated the identification of the relevant lipid species in the green alga *Chlamydomonas reinhardtii*. 106 species were identified within the lipid classes: phosphatidylserine, phosphatidylethanolamine, phosphatidylglycerol, phosphatidylinositol, monogalactosyldiacylglycerol, digalactosyldiacylglycerol, and sulfoquinovosyldiacylglycerol.

Novel Aspect

A tool for graphical analysis of HRMS/MS has been developed, using the Kendrick mass defect in combination with chromatographic characteristics.

References

1. Pluskal T., Castillo S., Villar-Briones A., Orešič M., BMC Bioinformatics, 11, 395 (2010).
2. Korf A., Vosse C., Schmid R., Helmer P. O., Jeck V., Hayen H., Rapid Communications in Mass Spectrometry, in press (2018).

ThP-116 / ANALYSIS OF POLYGLYCEROPHOSPHOLIPIDS USING ISOTOPE-LABELED METHYLATION BY NUPLC-ESI-MS/MS

Jong Cheol Lee - MyeongHee Moon

Yonsei University, Yonsei University/Chemistry, Seoul

Keywords: polyglycerophospholipids, (trimethylsilyl)diazomethane, nanoflow UPLC-ESI-MS/MS, Lipidomics, neurodegenerative disease

Introduction

Polyglycerophospholipids include phosphatidylglycerol (PG), bis(monoacylglycerol)phosphate (BMP), Hemi BMP, and bis(diacylglycerol)phosphate (BDP) having two glycerols, and cardiolipin (CL), monolysocardiolipin (MLCL), and dilysocardiolipin (DLCL) having three glycerols [1]. These polyglycerophospholipids are generally synthesized from PG as a precursor and closely related to each other. Therefore lipidomic analysis of these lipids is important.

Methods

Polyglycerophospholipids were methylated with (trimethylsilyl)diazomethane (TMSD) under an acidic condition. Deuterated reagents with DCI were used for relative quantitation by isotope-labeled methylation (ILM) method [2]. Efficacy of ILM for polyglycerophospholipids were evaluated using reversed-phase nanoflow UPLC-ESI-MS/MS.

Results

The molecular structures of methylated polyglycerophospholipid standards were investigated with the CID

experiments during nUPLC-ESI-MS/MS. The degree of methylation varied with the number of phosphate groups (1 for PG, BMP, Hemi BMP, and BDP, 2 for CL, MLCL, and DLCL). The efficiency of methylation for polyglycerophospholipids was found to be > 97% for most classes. Over-methylation in the free hydroxyl group of glycerol moieties were detected with less than 4.6%. PG and BMP which are geometrical isomers, and three isomers of BMP can be well resolved from each other after methylation by nUPLC-ESI-MS/MS. Furthermore, isomeric structures among Hemi BMP, MLCL, and DLCL can be well separated.

Conclusions

In this study, analysis of polyglycerophospholipids were investigated with ILM method using TMSD. As a result, ILM method with nUPLC-ESI-MS/MS enhanced the capability of analyzing polyglycerophospholipids qualitatively and quantitatively. The developed method can be utilized to analyze urinary polyglycerophospholipids from neurodegenerative diseases such as Parkinson's and Alzheimer's diseases.

Novel Aspect

Analysis of polyglycerophospholipids can be improved by using Isotope-labeled methylation (ILM) method in nanoflow UPLC-ESI-MS/MS.

References

1. Hullin-Matsuda, F., Luquain-Costaz, C., Bouvier, J., Delton-Vandenbroucke, I., Prostag. Leukotr. Ess., 81, 313-324 (2009).
2. Lee, J. C., Byeon, S. K., Moon, M. H., Anal. Chem., 89, 4969-4977 (2017).

ThP-117 / PLASMA LIPIDOMIC PROFILING OF FIVE DIFFERENT TYPES OF CANCER BY NUPLC-ESI-MS/MS

Gwnagbin Lee - Jong Cheol Lee - MyeongHee Moon

Yonsei Univ., Chemistry, Seoul

Keywords: cancer, biomarker candidate, lipids, nUPLC-ESI-MS/MS

Introduction

Lipids are the substances that play key roles in the metabolism of human bodies and are known to be related to cancer genesis pathways [1]. In this study, plasma lipids from 5 different cancer types (liver, stomach, lung, colorectal, and thyroid cancer) were investigated in order to discover potential lipid biomarkers of each cancer.

Methods

Lipids were extracted from patient's plasma using the modified Folch method with MTBE/methanol [2], and the extracted lipids were analyzed by nUPLC-ESI-MS/MS. Due to the study that blood lipid metabolites are affected by hormones in women after menopause, only male plasma samples were used in this study.

Results

As a result of analyzing the extracted lipids from patient's plasma with liver, stomach, lung, colorectal, and thyroid cancers along with that of healthy controls, a total of 357 lipids structures (58 PCs, 39 PEs, 26 PGs, 37 PIs, 26 PAs, 31 sphingo, 150 glycerol) were identified. Based on identified lipids, a total of 241 lipids (55 PCs, 58 PEs, 19 PGs, 24 PIs, 17 PAs, 21 sphingo, 47 glycerol) were quantified. Lipids showing significant changes (high abundance, >2-fold change, and one-way analysis of variance (ANOVA) with Bonferroni post hoc test with p-value<0.05) in each cancer type were selected: 31 liver, 16 stomach, 23 lung, 24 colorectal, and 16 thyroid cancer. Statistical analyses of these selected lipids were performed with logistic regression and receiver operating characteristic (ROC) curve and the area under the curve values (AUC) of each cancer were more than 0.8 except for thyroid cancer.

Conclusions

A total of 357 from patient's plasma with liver, stomach, lung, colorectal, thyroid cancer, and controls, were identified by nUPLC-ESI-MS/MS and among them 241 lipids were quantified by multiple reaction monitoring (MRM). Selected lipids showing significant differences from healthy controls can be utilized to develop characteristic biomarkers of cancers in the future.

Novel Aspect

Lipidomic profiles of five different types of cancer patient's plasma were systematically examined by nUPLC-ESI-MS/MS with MRM along with statistical evaluations.

References

1. Perrotti, F. Rosa, C. Cicalini, I. Sacchetta, P. Del Boccio, P. Genovesi, D. Pieragostino, D., *Int. J. Mol. Sci.*, 17, 1992-2017 (2016)
2. Bang, D. Y., Byeon, S. K., Moon, M. H., *J. Chromatogr. A*, 1331, 19-26 (2014)

ThP-118 / MALDI-TOF/MS ANALYSES DETERMINED URINARY PHOSPHOLIPIDS AS A POTENTIAL BIOMARKER FOR PROSTATE CANCER

Xin Li⁽¹⁾ - *Kenji Nakayama*^(1,2) - *Takayuki Goto*⁽¹⁾ - *Shusuke Akamatsu*⁽¹⁾ - *Koji Shimizu*⁽³⁾ - *Osamu Ogawa*⁽¹⁾ - *Takahiro Inoue*⁽¹⁾

Department of Urology, Kyoto University, Kyoto (1) - Support Center for Precision Medicine, Shimadzu Techno-Research Inc., Kyoto (2) - Clinical Research Center for Medical Equipment Development, Kyoto University Hospital, Kyoto (3)

Keywords: MALDI-TOF/MS, Urine; Prostate cancer, Phosphatidylcholine, Lysophosphatidylcholine

Introduction

The variations of phospholipids in plasma or tissues have been reported in several cancers, however they have not been elucidated in prostate cancer (PCa) especially in the non-invasive source of urine [1]. Matrix-assisted laser desorption and ionization time-of-flight mass spectrometry (MALDI-TOF/MS) is increasingly used in the lipid research which enables a rapid screening of phospholipids [2].

Methods

Urinary lipidomic screening was performed by MALDI-TOF/MS (AXIMA Performance, Shimadzu) in the urine samples collected from 20 PCa and 20 benign prostatic hyperplasia (BPH) patients after digital rectal examination (DRE). Lipids were identified by MALDI-TOF/MS2 analyses. The compositions of phosphatidylcholine (PC) and lysophosphatidylcholine (LPC) were compared between the PCa and BPH groups.

Results

In this study we have focused on the DRE urine. The samples contain many kinds of lipids secreted from prostatic microenvironment, which enable the detection of secreted prostate products as potential sources of PCa specific biomarkers. The lipids were extracted from the DRE urine by the Bligh-Dyer method[3]. The urinary phospholipids including PC, LPC, phosphatidylethanolamine (PE) and phosphatidylglycerol (PG) species were detected by MALDI-TOF/MS and identified by MALDI-TOF/MS2 in a reflectron positive mode. Among them, PC and LPC species were relatively abundant and with good reproducibility in the DRE urine samples. The peak intensities of the synthetic PC and LPC molecules spiked in the urine were detected by MALDI-TOF/MS and linearly correlated with their concentrations. The compositions of PC and LPC in the urine samples were significant different between the PCa and BPH groups.

Conclusions

The urinary phospholipids are detectable by MALDI-TOF/MS and identified by MALDI-TOF/MS². The relatively quantitative analyses of PC and LPC concentrations in the DRE urine samples are feasible and with good reproducibility. The MALDI-TOF/MS analyses determined the alteration of the compositions of urinary PC and LPC and it may be used as a simple and non-invasive biomarker for the PCa.

Novel Aspect

We revealed the variations of urinary PC and LPC detected by MALDI-TOF/MS in PCa which may offer a simple and non-invasive method for PCa diagnosis.

References

1. Bandu, R., Mok, H. J. & Kim, K. P. *Mass Spectrom. Rev.* 37, 107-138 (2018).
2. Fuchs, B., Süß, R. & Schiller, J. *Lipid Res.* 49, 450-475 (2010).
3. Bligh, E. G. & Dyer, W. J. *Can. J. Biochem. Physiol.* 37, 911-917 (1959).

ThP-119 / INVESTIGATION OF ONE NOVEL OXIDIZED 1-PALMITOYL-2-ARACHIDONOYL-SN-GLYCERO-3-PHOSPHOCHOLINE MOLECULE ON HUMAN AORTIC ENDOTHELIAL CELLS BY LC-MS/MS

Yu-Qing Chen - Hsin-Pai Tsou - Ying-Tzu Shieh - Mine-Yine Liu

National Changhua University of Education, Department of Chemistry, Changhua, Taiwan

Keywords: Phospholipids, oxidation, human aortic endothelial cell, LC-MS/MS

Introduction

Phospholipids are the major components of cell membranes. Among them, 1-palmitoyl-2-arachidonoyl-sn-glycero-3-phosphocholine (PAPC) is especially important since oxidized PAPC products are closely associated with atherosclerosis. PAPC is composed of a polyunsaturated fatty acyl chain at the Sn2 position, and is highly susceptible to oxidation. Many of the oxidized PAPC products are bioactive.

Methods

An LC-MS/MS method was developed to analyze one specific ox-PAPC molecule ($m/z=846.5$) in human aortic endothelial cells [1-5]. The optimal HPLC mobile phase compositions and flow rate were selected. For MS/MS, the optimal sheath gas pressure, capillary temperature, collision energy and injection time were chosen. Endothelial cells were treated with oxidized LDLs, extracted by liquid-phase and solid-phase extractions and then analyzed by LC-MS/MS.

Results

Elucidation of the chemical structure of $m/z=846.5$ was based on the MS² and MS³ spectra [1-5]. Quantitation of $m/z=846.5$ in HAECs was based on the reconstructed ion chromatogram (RIC) and selected multiple-reaction monitoring (SRM) of LC-MS/MS analysis. An internal standard was used to construct calibration curves for measuring the concentration of $m/z=846.5$ in endothelial cells under oxidative stress.

Conclusions

An LC-MS/MS method in conjunction with liquid-phase and solid-phase extractions has been developed to separate and identify a specific ox-PAPC molecule ($m/z=846.5$) in human aortic endothelial cells. The most possible structure of the ion is PAPC 5-epoxy_{8-E/D}-IsoP₁₄-epoxy which is elucidated from its MS² and MS³ spectra. This new ox-PAPC molecular species has never been reported by other researchers so far.

Novel Aspect

A novel ox-PAPC molecule ($m/z=846.5$) in human aortic endothelial cells has been found. Its chemical structure has been identified by LC-MS/MS analysis.

References

1. Berdeaux O., Juaneda P., Martine L., Cabaret S., Bretillon L., Acar N., J. Chromatogr. A., 1217, 7738 (2010).
2. Hui S. -P., Taguchi Y., Takeda S., Ohkawa F., Sakurai T., Yamaki S., Jin S., Fuda H., Kurosawa T., Chiba H., Anal. Bioanal. Chem., 403, 1831 (2012).
3. Watson A. D., Subbanagounder G., Welsbie D. S., Faull K. F., Navab M., Jung M. E., Fogelman A. M., Berliner J. A., J. Biol. Chem., 274, 24787 (1999).
4. Fam S. S., Morrow J. D., Curr. Med. Chem., 10, 1723 (2003).
5. O'Donnell V. B., Biochim. Biophys. Acta, 1811, 818 (2011).

ThP-121 / LIPID PROFILING OF NON-STEROIDAL ANTI-INFLAMMATORY DRUGS (NSAIDS) ON PROSTATE CANCER CELL LINES WITH TOF-SIMS AND DESI-MS

Danielle J. McDougall⁽¹⁾ - Emrys Jones⁽²⁾ - Claire Hart⁽³⁾ - Michael Brown⁽³⁾ - Adam McMahon⁽⁴⁾ - Nicholas Lockyer⁽¹⁾

University of Manchester, Division of Chemistry, Manchester (1) - Waters Corporation, Research and Development, Wilmslow (2) - University of Manchester, Division of Molecular and Clinical Cancer Sciences., Manchester (3) - University of Manchester, Division of Informatics, Imaging and Data Sciences, Manchester (4)

Keywords: Prostate Cancer, SIMS, DESI, NSAIDs

Introduction

Prostate cancer (PCa) is the prime cause of cancer related death in Western men. 1 It is not possible to distinguish indolent and aggressive disease, leading to overtreatment and associated complications. NSAIDs can potentially prevent growth of PCa cells, offering an alternative treatment approach. 2 The response of NSAIDs on PCa cell lines varying in clinical aggression was studied using MSI.

Methods

The lipidomic fingerprint was explored post NSAID treatment at a series of concentrations and time points. For DESI (Waters, Xevo-G2-SX), 95% methanol was used at a spatial resolution of 50 μm . For ToF-SIMS (Kore Technology Ltd), the samples were analysed using micro-focussed Au_3^+ clusters. Due to the hypoxic nature of the prostate, this was repeated under different levels of oxygenation.

Results

The cell lines used for this study were LNCaP, LNCaP-C42 and LNCaP-C42B, which respectively increase in aggression. Post drug treatment, the proliferation assays demonstrated that trifluorinated ibuprofen (TFIB), ibuprofen (IBP), and orlistat, lead to a decrease in cell growth above 1 mM of IBP concentrations and 10 μM of Orlistat. ToF-SIMS and DESI analysis of the untreated cell line samples highlighted a number of differences in the ratios of lipid classes, including more abundant cholesterol species in the less aggressive LNCaP cells with phospholipids generally being more abundant in the more aggressive cell lines. After drug treatment, multivariate analysis highlighted further differences. Under varying levels of hypoxia, the untreated samples separate with principle component analysis, suggesting future work identifying which more accurately represents the prostate in vitro.

Conclusions

The administration of NSAIDs leads to a decrease in PCa cell growth irrespective of their phenotype. Untreated PCa cell lines can be differentiated using ToF-SIMS and DESI with multivariate data analysis of the lipid fingerprint. NSAIDs alter these profiles suggesting lipid involvement in the drug pathway. Here we show the potential of MSI technologies to generate therapeutic mechanistic data.

Novel Aspect

Using imaging MS techniques to examine the lipid profile of PCa cells with differing metastatic properties post drug treatment.

References

1. Prostate Cancer (C61): 1993-2014. European Age-Standardised Incidence Rates per 100,000 Populations, Males, UK, 2014.
2. Andrews, J.; Djakiew, D.; Krygier, S.; Andrews, P. Superior Effectiveness of Ibuprofen Compared with Other NSAIDs for Reducing the Survival of Human Prostate Cancer Cells. *Cancer Chemother. Pharmacol.* 50, 277–284 (2002).

ThP-122 / DEFINING STRUCTURAL DIFFERENCES IN CARDIOLIPINS FROM AN ACTINOMYCETE MARINE SPONGE SYMBIONT AND S. AUREUS

Benjamin Oyler⁽¹⁾ - ***Courtney Chandler***⁽²⁾ - ***Fan Zhang***⁽³⁾ - ***Robert Ernst***⁽²⁾ - ***Russell Hill***⁽³⁾ - ***David Goodlett***⁽⁴⁾

University of Maryland, Baltimore, School of Medicine, Baltimore (1) - University of Maryland, Baltimore, School of Dentistry, Baltimore (2) - University of Maryland Center for Environmental Science, Institute of Marine and Environmental Technology, Baltimore (3) - University of Maryland, Baltimore, School of Pharmacy, Baltimore (4)

Keywords: cardiolipin, mass spectrometry, ion mobility, symbiont, actinomycete

Introduction

Marine sponges harbor complex microbiomes and novel bacteria associated with these sponges have been isolated and classified by 16S rDNA sequencing [1]. Sponge-symbiont relationships are not yet well understood. It is clear, however, that symbionts are responsible for a portion of the production of bioactive metabolites previously isolated from sponges. Some membrane lipids also have been shown to play a vital role in bacterial recognition by host cells.

Methods

An actinobacterium (R45601) was isolated from X. mutant cultured in ISP2 medium. Cells were lysed and cardiolipins (CLs) were extracted according to the method described by El Hamidi, et al [2] for extraction of LPS. Extracts were dissolved in a 2:1:0.1 solution of chloroform, methanol, and water. Samples were analyzed first by MALDI-TOF using norharmane as matrix. Tandem mass spectra with ion mobility separation were acquired on a Q-IMS-TOF mass spectrometer.

Results

A comparison of broadband mass spectra of extracted lipids, filtered by drift time in the ion mobility cell and m/z , from a clinical isolate [3] of methicillin-resistant *S. aureus* (MRSA-M2) and a recently discovered [4], marine sponge-associated actinomycete (R45601) showed that R45601 produces a pool of CLs with greater mean mass and possibly a higher degree of unsaturation than those produced by MRSA-M2. To identify fatty acyl groups incorporated into CLs from each species, MS/MS experiments were performed. The fatty acyl groups from CLs were released in the gas phase upon collisional activation as ionized free fatty acids with a single negative charge. These data show that ions observed in the single stage mass spectra which differ by multiples of 2 Da correspond to CLs with unsaturated fatty acyl groups. CLs from MRSA-M2 contained mostly saturated fatty acyl groups, whereas many CLs from R45601 contained at least one unsaturation. A possible explanation for this difference is that the mechanisms for CL synthesis differ between bacterial phyla.

Conclusions

Comparisons of molecular classes from different organisms can be useful for generation of testable

hypotheses. In this study, there were clear observations of differences in CL structure between representatives from two bacterial phyla. These structural differences may have implications with regard to bacteria-host interactions. CLs may play a role in sponge recognition of microbes. Differences in CL structure will also affect bacterial phenotype and behavior.

Novel Aspect

Direct comparison of cardiolipins from a terrestrial firmicute and a marine actinomycete.

References

1. Hentschel U., Piel J., Degnan S., Taylor M. Nat. Rev. Microbio., 10, 641 (2012).
2. Hamidi E., Tirsoaga A., Novikov A., Hussein A., Caroff M. J. Lipid Res. 46, 1773 (2005).
3. Shirliff M., Calhoun J., Mader J. Clin. Orthop. Relat. Res. 401, 239 (2002).
4. Montalvo N., Davis J., Vicente J., Pittiglio R., Ravel J., Hill R. PLoS One. 9, e90517 (2014).

ThP-123 / UNTARGETED LIPIDOMICS FOR THE ASSESSMENT OF ARSENIC EXPOSURE IN RICE (ORYZA SATIVA JAPONICA)

Miriam Carolina Pérez Cova - Romà Tauler Ferre - Joaquim Jaumot Soler

IDAEA-CSIC, Department of Environmental Chemistry, Barcelona, Spain

Keywords: untargeted, lipidomics, LC-MS, arsenic, rice.

Introduction

There is an increasing interest in how arsenic is assimilated and metabolized in plants-base foods such as rice (*Oryza sativa japonica*), due to the potential health risk that its accumulation might cause through food chain[1]. Lipidomics is a powerful approach to study lipids species present in rice, with the aim of determining the affected metabolic pathways.

Methods

Two treatments were applied to rice crops: As was applied by irrigation (1 or 1000 μM), and directly to the soil before planting(5 or 50ppm). LC-MS analysis was conducted using a Waters Acquity UPLC system connected to a Waters TOF mass spectrometer. A Kinetex C8 column (100 x 2.1 mm, 1.7 μm) was used for chromatographic separation. Multivariate tools (PCA, PLS-DA, ASCA) were employed in data analysis.

Results

In watering treatment, the lower concentration was set at 1 μM , as it is the limit of acceptable As concentration in water by European legislation, whereas the higher concentration was chosen at 1000 μM , because it is expected to induce changes in plants without causing their death [2]. In soil treatment, the lower As limit was set at 5 ppm, as an average concentration of this element in Earth's crust, and the higher As limit at 50 ppm, as soil is considerate to be contaminated by this element above 10 ppm [1].

Results obtained after chemometric analysis show that there are significant differences between treatments and concentrations, when compared with control samples (without arsenic treatment), especially in the case of the higher concentrations of this element.

The principal lipid families that were affected by As treatments were glycerophospholipids and acylglycerides.

Conclusions

Only the highest As concentration levels (both in irrigation water or settled in soil) allowed the differentiation between samples. This complies with the European legislation, with permitted levels of the As exposure not showing relevant effects.

As a general lipidomic conclusion, glycerophospholipids and acylglycerides are the main families of lipids that were affected by the arsenic exposure.

Novel Aspect

A lipidomic approach giving a better understanding of the metabolic pathways affected by arsenic exposure, and able to identify possible biomarkers.

References

1. Zhao FJ., McGrath SP., Meharg AA., *Annu. Rev. Plant Biol.* 61, 535-59 (2010).
2. Navarro-Reig M., Jaumot, J., García-Reiriz, A., Tauler R., *Anal. Bioanal. Chem.* 407, 8835-8847 (2010).

ThP-124 / HIGH-RESOLUTION MULTIPLE REACTION MONITORING METHOD DEVELOPMENT FOR QUANTIFICATION OF STEROIDAL HORMONES IN PLASMA

Ana Paula Ferranti Peti⁽¹⁾ - ***Gisele Aparecida Locachevic***⁽¹⁾ - ***Morgana Kelly Borges Prado***⁽¹⁾ - ***Luiz Alberto Beraldo de Moraes***⁽²⁾ - ***Lúcia Helena Faccioli***⁽¹⁾

Universidade de São Paulo, FCFRP, DACTB, Ribeirão Preto (1) - Universidade de São Paulo, FFCLRP, Chemistry Department, Ribeirão Preto (2)

Keywords: Lipidomics; MRM-HR; LC-MS/MS; Steroids; Chronic Stress

Introduction

Advances in MS enabled development of HR-MRM and became the analysis of target analytes more specific in lipidomics approach, contributing in the medical and pharmaceutical fields. Hence, HR-MRM it is a sensitive and very specific virtual modes of acquisition that allows the simultaneous steroidal hormones quantification in low concentrations, such as observed in biological samples.

Methods

HR-MRM method included cortisol, cortisone, corticosterone, 11-DHC, progesterone and aldosterone. Analytical validation was conducted according to the ANVISA guidelines. Specificity, linearity, intraday and interday precision and accuracy, and the LLOQ were evaluated. Plasma real samples of 129 sv male mice exposed to chronic unpredictable stress were used to validate the quantification method.

Results

The method was conducted according to the ANVISA normative, adopting a coefficient of variation (CV), as well as relative standard deviation (RSD) and relative error (RE) lower than 15% in linearity, intraday and interday precision and accuracy. For cortisol, corticosterone, and their inert metabolites (cortisone and 11-DHC), the LLOQ was 3.9 ng/mL, while those for progesterone and aldosterone was 7.8 and 15.6 ng/mL, respectively. The linearity of the method was between LLOQ-500 ng/mL for each analyte. HR-MRM analysis showed that animals submitted to stressors have 4.5 times more corticosterone in their serum than non-stressed mice and the levels of 11-DHC does not vary significantly in both groups. The development method showed adequate selectivity and sensibility for quantification of target analytes and can be used in several studies according the observed validation parameters.

Conclusions

The results indicate that the method can be applied for quantification of steroids in several biological samples, such as human or animal plasma with selectivity and confidence. Thereby all studies that involve the balance of the target analytes can be conducted using this method for quantification.

Novel Aspect

The HR-MRM method allows the monitoring of different steroids from human/animal origin, facilitating its

application in several biological situations.

References

1. Peti A. P. F., Locachevic G. A., Prado M. K. B., Moraes L. A. B., Faccioli L. H., J. Mass Spectrom., 53, 423-431 (2018).
2. Stewart P. M., Mason J. I., Steroids, 60, 143-146 (1995).
3. Galvão A. F., Petta T., Flamand N., Bollela V. R., Silva C. L., Jarduli L. R., Malmegrim K. C., Simões B. P., Moraes L. A., Faccioli L. H., Anal. Bioanal. Chem., 408, 3613-3623 (2016).

ThP-125 / MICRO-UPLC-MS HIGH-THROUGHPUT SCREENING ASSAY FOR SPHINGOLIPIDS PATHWAY ANALYSIS

Kristen Randall⁽¹⁾ - ***Helen Klodnitsky***⁽¹⁾ - ***Drew Tietz***⁽¹⁾ - ***Clifford Phaneuf***⁽¹⁾ - ***Maureen Olszewski***⁽²⁾ - ***Elina Makino***⁽²⁾ - ***Keeley Murphy***⁽³⁾ - ***Alexei Belenky***⁽¹⁾

Sanofi, Pre-Development Sciences, Waltham (1) - Sanofi, In Vitro Biology, Waltham (2) - Thermo Fisher, Marketing, San Jose (3)

Keywords: microflow liquid chromatography, high-throughput screening, lipidomics

Introduction

In the presentation we will describe the Micro-UPLC-MS system based on Vanquish Focused+ LC (Thermo) and TSQ Altis Triple Quadrupole Mass spectrometer (Thermo) for lipid analysis.

Methods

Replacing conventional tubing of the system with micro ID tubing, adding a second pump and optimizing the injection strategy allows us to convert conventional system to micro with parallel chromatography capabilities. Run time using a 5 cm length column was less than one min per injection.

Results

This system was successfully used for affinity screening of biomarkers from lipid pathways. This strategy is valuable when highly sensitive biomarker quantification in large sample set (library screening) is required.

Conclusions

In the poster we will demonstrate performance of the micro-LC-MS system using endogenous sphingolipid biomarkers from cell extracts at sub ng/ml (limit of quantitation is 0.05 ng/ml) concentration. We have evaluated Acclaim C8 0.3 mm i.d. and BEH C18 1 mm i.d. columns using gradient and isocratic chromatography. C8 was found to be more rugged but C18 enabled us to discriminate isobaric lipids.

Novel Aspect

We have developed LC-MS based screens using microflow chromatography.

ThP-126 / CERAMIDOME PLASTICITY THROUGH ISO-ENERGETIC PRECURSOR AND NEUTRAL LOSS DISCOVERY SCANS

Federico Maria Rubino - ***Michele Dei Cas*** - ***Jessica Rizzo*** - ***Paola Signorelli*** - ***Riccardo Ghidoni*** - ***Rita Clara Paroni***

Universita degli Studi di Milano, Dipartimento di Scienze della Salute, Milano

Keywords: Ceramide, epidermoside, collision-induced dissociation, screening, triple quadrupole,

Introduction

Mass spectrometric methods of identification and measurement are pivotal to investigate the multi-faceted biological role of ceramides and the dose-effect relationship of pharmaceuticals that modulate ceramide flux. Triple-quadrupole tandem MS embodies the highest flexibility to discover, identify and measure expected and unexpected molecular species, and an unprecedented scan mode is presented [1].

Methods

An API3000 LTQ LC-MS system is operated under manufacturer's instructions to perform extensive tandem MS fragmentation studies on standard ceramides with C₆-C₂₄ fatty acids, lipid extracts of animal tissues and of nutraceutical seeds and nuts. LC separation and MS operation is optimized for the separation, detection, identification and quantification of expected and unexpected ceramides.

Results

The equation linking instrumental to center-of-mass collision energy (CoM-CE) in collision-induced dissociation is rearranged to calculate laboratory-frame inter-quadrupole potential drop (LF-CE) corresponding to distinctive positions, such as that of the CoM-CE maximum fragment yield, in the breakdown curves for protonated ceramides with molecular masses up to 1,000 Th. The value of the CoM-CE for distinctive fragments is calculated for representative ceramide FA homologs from accurate recordings of the spectra in the 0 to 90 DV interval of LF-CE. A scan function of the q2 collision cell is synchronized with the Q1 scan in the precursor ion mode and with the Q1 and Q3 scans in the Neutral Loss mode. Several lipid extracts containing ceramides are examined by infusion and by LC, and unexpected components, such as Z-ceramides in nut extracts are tentatively identified.

Conclusions

Precursor and Neutral Loss scans with iso-energetic CID in the triple quadrupole MS allow to characterize ceramide mixtures that contain unexpected components with a common substructure to known ones. Individual response factors for (semi)-quantification of unexpected components can be calculated from the fractional intensity of the reporter fragment in the spectra of known components.

Novel Aspect

A novel scan function of a triple quadrupole instrument allows to measure tandem mass spectra at a constant center-of-mass collision energy.

References

1. Rubino F. M., Campisi G. M., Paroni R., Proc. XII Sphingolipid Club Meeting, Trabia (Italy), Sep. 6-10, 2017 T7, P6 (2017).

ThP-127 / SERUM LIPIDOMIC CHANGES IN ZIKA AND DENGUE VIRUS INFECTION REVEALS POTENTIAL TARGETS FOR ASSERTIVE CLINICAL DIAGNOSIS

Fabio Santos⁽¹⁾ - ***Aline Martins***⁽²⁾ - ***Kelly Magalhães***⁽²⁾ - ***Marcos N. Eberlin***⁽³⁾

University of Campinas, Analytical Chemistry, Campinas (1) - University of Brasilia, Cell Biology, Brasilia (2) - University of Campinas, Organic Chemistry, Campinas (3)

Keywords: Zika virus, Dengue virus, lipidomic, serum

Introduction

The Zika (ZKV) and Dengue (DENV) Virus have been a major public health challenge in the world due to

large number of underreporting, inaccurate and false diagnoses. The lipidomic-based method to detect biomarkers of the physiological variations at all stages of the infection and immunological responses. Herein we showed that lipidomics-based mass spectrometry can be applied for molecular diagnostic of Zika and Dengue and co-infection of infected serum [1,2].

Methods

The serum samples were separated in four groups through ELISA (IgG and IgM) protocol: 30 ZKV, 30 DENV and 30 controls. The method of metabolites extraction from 50 μ L serum infected was Bligh-Dyer. The upper phase was collected for lipidomic analysis using an Agilent 1290 Infinit LC-MS system coupled with an Agilent 6550 Q-TOF/MS mass spectrometer equipped with an electrospray source, in positive and negative modes.

Results

The chromatograms were processed by XCMS software and multivariate analysis was made in the Metaboanalyst 3.0 software. The PLS-DA was applied to classify the serum samples according to the type of virus infection. Based in different class of lipids, the serum infected by ZKV and DENV were clearly distinguishable into two clusters by PSL-DA plot. The majors lipids identified as responsible for the differentiation of the two clusters, by variable importance in projection (VIP plot), were phosphatidylcholine (PC), lysophosphatidylcholine (LysoPC), phosphatidylinositols (PI), glycolipids and glycosphingolipids. The Gangliosides, particularly GM1, GM2, GD1a, GD1b, and GT1 were upregulated only under infection of the ZKV.

Conclusions

Our study identified altered lipids phosphatidylcholine (PC), lysophosphatidylcholine (LysoPC), phosphatidylinositols (PI), glycolipids and glycosphingolipids of Zika and Dengue infection. The identification of specific biomarkers of the Zika and Dengue virus infections can perform the differential diagnosis of these arboviruses and closer to the actual clinical state of the patient.

Novel Aspect

The lipids identified allude to possible novel diagnostic targets, pharmacological studies and advance our understanding on the mechanisms of arbovirus pathogenesis.

References

1. Melo C. F. O. R., Delafiori J., de Oliveira D. N., Guerreiro T. M., Esteves C. Z., Lima E. D. O., Pando-Robles V., Catharino R. R. *Front. Microbiol.* 8, 1954 (2017).
2. Cui L., Pang J., Lee Y. H., Ooi E. E., Ong C. N., Leo Y. S., et al. *PLoS Negl Trop Dis* 12(1): e0006217 (2018).

ThP-128 / FISHING FOR LACTONES: METHOD DEVELOPMENT FOR SELECTIVE IDENTIFICATION OF LACTONE-CONTAINING LIPIDS

Elana Slutsky Smith⁽¹⁾ - *Soliman Khatib*⁽²⁾ - *Andrea Szuchman-Sapir*⁽¹⁾

Migal-Galilee Research Institute, Laboratory of Vascular Signaling Research, Kiryat Shmona (1) - Migal-Galilee Research Institute, Oxidative Stress and Human Diseases Laboratory, Kiryat Shmona (2)

Keywords: LC-MS, lactones, PON1, polyunsaturated fatty acids (PUFA), structure determination

Introduction

PON1 is an anti-atherogenic HDL-associated lactonase, whose native substrates remain unknown. PON1 knockout mice have low blood pressure and high levels of arachidonic acid (AA) metabolite 5,6-EET [1]. Under physiological conditions, 5,6-EET transforms to the δ -lactone 5,6- δ -DHTL [2], a PON1 substrate [3].

We hypothesize that additional lactone PUFA metabolites are potential PON1 substrates.

Methods

To identify and characterize these lactones, we react lipid extracts from mouse kidneys with an amine to convert lactones to lactams that can be isolated using solid phase cation exchange. We apply high resolution liquid chromatography-mass spectrometry to determine the molecular formula of the original lactones. Molecular databases are used to find potential structures, and standards will be used to confirm or negate these structures.

Results

The selective chemical reaction was developed using standards of different lipid functional groups (free fatty acid, ester, and two lactones of different ring size). Conditions were optimized to achieve reaction of the lactones and not of the ester. However, reaction of the fatty acid did occur. LC-MS analysis of similarly reacted lipid extracts from w. t. and PON1KO mice provided over 60 different chemical formulae of potential lactones, most of which were present in all samples (n = 3 pools of 5 kidneys for each strain). Using molecular databases, hundreds of potential lactone structures were identified. We are currently investigating ways to limit the search results to the most biologically relevant and likely structures. Standards of those structures will then be subjected to the same chemical reaction, extraction, and LC-MS method in order to confirm the structures of the unknowns in the biological samples. MS/MS experiments may be applied. It will then be investigated whether the confirmed structures are PON1 substrates.

Conclusions

Chemical formulae of potential lactones were repeatedly identified in separate mouse kidney samples using LC-MS, while their precise structures are not yet known. Presently, further repetitions are underway to determine the differences between w. t. and PON1KO and pinpoint PON1 substrates. We hope that this research will shed light on the potential role of PON1 in blood pressure regulation.

Novel Aspect

PON1 has been thoroughly studied in relation to atherosclerosis, yet the substrates of PON1 linking the enzyme to blood pressure regulation remain unknown.

References

1. Gamliel-Lazarovich A., Abassi Z., Khatib S., Tavori H., Vaya J., Aviram M., Keidar S., *Atherosclerosis*, 222, 92 (2012).
2. Fulton D, Falck J. R., McGiff J. C., Carroll M. A., Quilley J., *J. Lipid Res.*, 39, 1713 (1998).
3. Eryanni-Levin S., Khatib S., Levi-Rosenszvig R., Tamir S., Szuchman-Sapir A., *Biochim. Biophys. Acta*, 1851, 1118 (2015).

ThP-129 / CHALLENGES OF LIPID ISOMERS AND ISOBARS IN MASS SPECTROMETRY-BASED BIOMARKER RESEARCH

Ulf Sommer - Hai Pham-Tuan - Therese Koal

Biocrates Life Sciences AG, Innsbruck, Austria

Keywords: Isobars and isomers, lipidomics, high-resolution accurate-mass (HRAM) spectrometry, multiple-reaction monitoring (MRM)

Introduction

Interfering isobaric & isomeric species still constitute a major issue in the unambiguous attribution of lipid biomarkers to disease. MRM scanning allows to distinguish the main lipid classes, but not always isobaric subclasses. HRAM spectrometry (MS1) allows the distinction of these isobaric species but is not able to

distinguish isomers. We are comparing two flow-injection analysis methods.

Methods

Two commercially available kits using MRM or HRAM screening are being compared. They are used to measure a range of small molecules & lipids and have been validated for human plasma & serum, the AbsoluteIDQ® p180 Kit (1-3) for Waters Xevo™TQ-S, Sciex QTrap®5500 or Thermo Scientific TSQ Vantage™, and p400 HR Kit (2017)(ThermoQ Exactive™ systems) (both Biocrates Life Sciences AG).

Results

Derivatization removes especially most phosphatidylethanolamines, the main interfering species, leading to more specific detection of PCs.

All realistically possible species covered by a result should be listed. The given name of the signal should refer to the most relevant species; lists with isobaric / isomeric species being provided.

In both cases, chromatography will assist the separation of species, but is limited due to high throughput requirements and the need for internal standards for different sections of the chromatogram to correlate with suppression effects.

The reliable quantitation of additional lipid classes can be achieved, *e.g.* as a service (4), but requires additional effort compared to the single high-throughput methods discussed here.

Conclusions

Screening using HRAM and low-resolution MRM will by their nature give somewhat different results due to different isobaric & isomeric compounds. The AbsoluteIDQ® p180 and p400 HR Kits provide standardized procedures for the fast and efficient detection and annotation of the major lipid classes. Use of HRAM increases coverage through complementary quantitative high-resolution measurements.

Novel Aspect

New guidance on how to interpret the different results obtained by lipid biomarkers screening via MRM versus HRAM.

References

1. Draisma H. H. M. et al., Nat. Commun. 6, 7208 (2015).
2. Klavins K. et al., Alzheimers Dement. 1, 295-302 (2015).
3. Siskos A. P. et al., Anal. Chem. 89, 656-665 (2017).
4. Köberlin M. S. et al., Cell 162, 170-183 (2015).

ThP-130 / APPLYING TRAPPED ION MOBILITY SEPARATION (TIMS) IN COMBINATION WITH PARALLEL ACCUMULATION SERIAL FRAGMENTATION (PASEF) FOR ANALYSIS OF LIPIDOMICS SAMPLES

Matthias Szesny⁽¹⁾ - ***Sebastian Götz***⁽¹⁾ - ***Sven Meyer***⁽¹⁾ - ***Ulrike Schweiger Hufnagel***⁽¹⁾ - ***Aiko Barsch***⁽¹⁾ - ***Ningombam Sanjib Meitei***⁽²⁾

Bruker Daltonik GmbH, Metabolomics/Lipidomics, Bremen (1) - Premier Biosoft, Lipidomics, Palo Alto (2)

Keywords: PASEF; Ion mobility; Lipidomics; timsTOF PRO; Plasma

Introduction

Untargeted lipidomics workflows are aimed at profiling changes in the lipidome which can lead to the discovery of lipids as potential biomarkers. We present an improved lipidomics identification workflow that is based on the combination of LC-MS and TIMS separation together with a data dependent MS/MS fragmentation (PASEF).

Methods

LC separation was performed using a Bruker Elute UHPLC system. Data acquisition was performed with separate positive and negative ionization methods on a Bruker timsTOF PRO instrument. The instrument was in TIMS mode while acquiring MS/MS fragmentation (PASEF). The data sets were processed in MetaboScape 4 using 4-dimensional feature finder. Features were annotated using SimLipid software.

Results

The TIMS technology adds an additional dimension of separation. To evaluate this effect on lipidomics analyses a direct infusion method was applied to critical pairs of isobaric and close isobaric lipids. A pair of isobaric phosphocholines (PC 18:1) differing only in the position of the unsaturation (9Z vs 6Z) could be separated with a mobility peak resolution of 175. A second almost isobaric pair (phosphoethanolamine & phosphocholine) could be fully baseline separated with a mobility resolution of 180.

LC-TIMS-MS measurements in PASEF mode were performed on extracts of commercially available plasma samples from different species (bovine, chicken, pig and human). The 4-dimensional feature finder in MetaboScape was used to extract the features from all analyses. The main precursor mass and the most intense fragment spectrum were exported to SimLipid software for identification.

Conclusions

First results show 200-400 unique lipid IDs (depending on animal species). Compared to analyses performed without TIMS/PASEF an increase of 30-50% in identifications can be observed. Additional measurements will be presented to further highlight the benefit of PASEF.

Novel Aspect

Lipidomics workflow combining positive/negative ion data using PASEF mode on an LC-TIMS-QTOF instrument.

ThP-131 / HOW TO ENSURE QUALITY IN THE ANALYSIS OF ENDOGENOUS COMPOUNDS IN BIOLOGICAL SAMPLES

Dominique Thomas⁽¹⁾ - ***Stefanie C. Fleck***⁽²⁾ - ***Robert Gurke***⁽¹⁾ - ***Gerd Geisslinger***⁽¹⁾ - ***Nerea Ferreirós***⁽¹⁾

Goethe University, Institute of Clinical Pharmacology, Frankfurt (1) - Fraunhofer IME, Project Group TMP, Frankfurt (2)

Keywords: endogenous compounds, quality assurance, validity of results

Introduction

The mass spectrometric analysis of lipid mediators and other endogenous compounds remains challenging and rich in pitfalls. Variations between runs and even within runs impairing the validity of the data are hard to identify. We established a systematic approach to ensure the quality of data in the analysis of prostanoids in plasma samples by evaluating control samples using a control chart.

Methods

A mass spectrometric method for the quantification of PGE₂, PGD₂, TXB₂, PGF₂α and 6-keto-PGF₁α in human plasma was validated according to EMA and FDA guidelines. After successful validation, a control chart was set up: blank plasma was spiked with prostanoids and aliquots were stored at -80°C. In every analytical run, one control plasma aliquot was analysed in the beginning and one at the end of the batch, in addition to six freshly prepared quality control samples.

Results

After analysis of 24 control samples over six weeks, for all analytes alert limits and control limits were

defined depending on mean value and corresponding standard deviation. Exemplarily, results for PGE2 were: MV 93.2 pg/mL, SD 8.7 pg/mL, alert limit $\pm 2 \cdot SD$ (75.8 pg/mL or 110.6 pg/mL), control limit $\pm 3 \cdot SD$ (67.1 pg/mL or 119.3 pg/mL). After six more weeks, a tendency of increased concentrations of all prostanoids in the control samples was observed, indicating instability of the analytes in plasma samples at -80°C after 12 weeks. During validation, stability had been proven for 8 weeks.

Surprisingly, concentrations went back to the initial level after 10 more weeks. Detailed investigation revealed that the slopes of the calibration curves of the runs with increased concentrations of control samples differed from the slopes during the validation. These variations cannot be discovered when using only quality control samples. Thus, the control chart revealed difficulties in the preparation of calibration and quality control standard solutions.

Conclusions

Due to lack of interlaboratory surveys for lipid mediators and other endogenous compounds, alternative tools for quality assurance are urgently required. Monitoring of control samples using control charts is a useful tool to ensure the constancy of the analytical results. Careful evaluation can reveal systematic errors, problems concerning stability of samples or analytical standards, lack of robustness and other issues relevant to quality.

Novel Aspect

A systematic approach to ensure quality and validity in the analysis of endogenous compounds.

ThP-132 / STUDY OF ACETOGENINS IN A NATURAL EXTRACT USING SFC-HRMS/MS AND POST-COLUMN METAL CATIONISATION

Laurent Laboureur⁽¹⁾ - *Pierre Champy*⁽²⁾ - *David Touboul*⁽¹⁾

CNRS, ICSN, Gif-sur-Yvette (1) - CNRS UMR 8076 BioCIS, Pharmacy Faculty, Chatenay-Malabry (2)

Keywords: supercritical fluid chromatography, high resolution tandem mass spectrometry, post-column metal cationisation, lipid.

Introduction

Acetogenins, lipids from the polyketide class and specifically found in plants from the Annonaceae family, are food toxins suspected to be responsible of atypical Parkinsonian syndromes. Acting as an inhibitor of mitochondrial complex I, acetogenins are also proposed as potential anti-cancer molecules. Previous dereplication studies highlighted some limitations of the HPLC-MS approach including long analysis time and low chromatographic resolution.

Methods

Nowadays, supercritical fluid chromatography (SFC) can be considered as an alternative to LC for the study of hydrophobic compounds. Indeed, physicochemical properties of supercritical CO_2 allow excellent solubility of such analytes and give analytical performances above that of LC. However, examples using SFC-HRMS/MS to analyze natural extracts are poorly developed in the literature.

Results

Due to the rarity of standards, a SFC-HRMS/MS method was directly developed using a crude seed extract of *Annona muricata* L. with a reduced analysis time by a factor of 4. Nevertheless, the fragmentation of protonated or sodium/potassium cationized species only led to limited structural information. For the first time, post-column lithium cationisation was fully optimized after SFC in order to form ion fragments of structural interest but of very low intensity. Finally, copper salt was employed allowing the detection of intense odd and even ion fragments related to unique radical-driven fragmentation pathways.

Conclusions

The final SFC-HRMS/MS method significantly improved the structural characterization of acetogenins highlighting the power of SFC-HRMS/MS in the field of natural products.

Novel Aspect

For the first time, SFC-HRMS/MS including post-column cationization was implemented for the structural characterization of natural lipids.

ThP-133 / IMPROVED LIPID ANNOTATION UTILIZING POSITIVE AND NEGATIVE ION MS2 / MS3 HCD AND CID SPECTRA

David A. Peake⁽¹⁾ - *Reiko Kiyonami*⁽¹⁾ - *Daniel Gachotte*⁽²⁾ - *Yasuto Yokoi*⁽³⁾ - *Sally Webb*⁽¹⁾ - *Andreas Huhmer*⁽¹⁾

Thermo Fisher Scientific, CMD Life Sciences Research, San Jose (1) - Dow AgroSciences, Indianapolis (2) - Mitsui Knowledge Industry, Tokyo (3)

Keywords: high-resolution, accurate-mass, annotation, LC-MSn, LipidSearch

Introduction

We present improvements to Thermo Scientific LipidSearch software. New algorithms are introduced to reduce false positives, improve quantitation and to automate searching of LC-MSn data obtained by higher collisional energy (HCD) and linear ion trap collisional induced dissociation (CID) fragmentation methods. The use of LC-MSn is applied to the more complete characterization of triacylglycerols in total lipid extracts from bovine tissue and insect larvae.

Methods

Bovine liver and insect larvae total lipid extracts were analyzed by data dependent LC-MS2 and MS3 using a Tribrid Orbitrap instrument. The MSn data were searched against the m/z of selected lipid precursor ions in the LipidSearch database and their predicted product ions and neutral losses. Each lipid adduct ion is ranked by mass tolerance, match to the predicted fragmentation and fraction of total MS-MS intensity.

Results

The number of lipid species annotated in each experiment were assessed at the sum composition (MS) and isomer (MSn) levels. LC-dd-MS2/MS3 spectra for potential lipid species were annotated separately from positive and negative ion adducts. Each lipid annotation was correlated within a chromatographic time window by merging the positive and negative ion MSn data into a single lipid result. This approach provides lipid annotation that reflects the appropriate level of MS2/MS3 product ions and neutral losses from the entire dataset giving higher confidence in lipid annotations. The merged results were filtered by minimum number of data points, signal-to-noise ratio, adduct ion, match score, ID quality, and coefficient of variation from replicate sample injections. Compared to the results generated only from dd-MS2 HCD results, a combination of HCD and CID LC-MSn gave significantly higher quality lipid identifications. From bovine liver the number of PC and TG results with complete acyl chain information improved by 12% to 23%, respectively.

Conclusions

Lipid annotation is often over-reported because the current software/database approaches match mass spectral evidence to exact lipid structures [1]. The lipidomics community is working towards reporting the correct level of annotation based upon the available mass spectral information. Utilizing a combination of HCD MS2 and intelligent product ion or neutral loss triggered CID MS2/MS3 data improves lipid annotation available within a single LC/MSn experiment.

Novel Aspect

LipidSearch 4.2 new software for improved lipid annotation and increased confidence using positive/negative ion, HCD/CID LC/MSn data from Orbitrap mass spectrometers

References

1. RustamY. H., Reid, G. E. *Anal. Chem.*, 90(1), 374-397 (2018).

ThP-134 / A SIMPLE AND RAPID EXTRACTION OF LIPIDS IN PLASMA USING SPIN COLUMN WITH SUPERABSORBENT POLYMER BEADS

Geul Bang⁽¹⁾ - *Jung Hoon Choi*⁽¹⁾ - *Young Ja Hwang Park*⁽²⁾ - *Jeong Ah Kim*⁽¹⁾ - *Young Hwan Kim*⁽¹⁾

Korea Basic Science Institute, Biomedical Omics Group, Cheong-ju, Republic of Korea (1) - Korea University, College of Pharmacy, Sejong, Republic of Korea (2)

Keywords: lipid extraction, superabsorbent polymers, spin column, liquid chromatography, mass spectrometry

Introduction

A fast and straightforward lipid extraction method from plasma is introduced in this study. Lipids are one of the major components of biological cell membranes, and they participate in many biochemical functions and metabolic processes, such as intercellular signaling, secretion, and energy storage [1]. The aim of our study was to develop a lipid extraction using a spin column with superabsorbent polymer beads (mSAPs). Although our previous study about fast lipid extraction method using SAPs has been reported [2], it is still challenging to realize quantitative analysis. Thus, a newly developed a simple and reliable lipid extraction method using a spin column filled with SAPs was built to shorten the time of extraction. We evaluated the extraction efficiency and the reproducibility of this method using liquid chromatography-mass spectrometry (LC/MS) analysis to compare properties of our method with the modified Folch method. The modified SAP method was applied to lipid extracts from the plasma sample, demonstrating that it can be powerfully utilized for high-speed (<10 min) preparation for lipid extraction.

Methods

The ten microliter of plasma was dropped into the spin column with SAP beads and left for 30 s to complete absorption of aqueous solutions. Then, mixture of organic solvents, MTBE:MeOH, 2:1 (v/v) were loaded to the SAP beads, and the samples were incubated for 2–3 min. For obtaining extracts, the spin columns are followed by centrifugation at 6,000 rpm for 2 min. The organic phase was collected and dried with nitrogen gas. The extracted lipids were dissolved in 200 μ L of IPA:ACN:DW (65:30:5 v/v) for storage at 4°C before MS analysis. Then, they were analyzed using ultra-performance liquid chromatography–mass spectrometry (UPLC/MS).

Results

There is an urgent requirement of a quantitative lipid extraction method for plasma sample especially in the area of clinical research. For the conventional lipid extraction, arduous supernatant transfer and large volume organic solvents need to be repeated twice. Hence, we have to adopt a spin column with SAPs. After its improvement, time-consuming reduces about 95 % for lipid extraction. For the qualitative lipid profiling after extraction using the modified SAPs method, plasma samples were analyzed using LC-ESI/MS. The results of the LC/MS chromatograms and MS/MS spectra indicated that the modified SAPs method successfully extracted lipids from plasma. Also, we found that LOD of the modified SAPs method (0.035pmol) was approximately about seven times lower than that of the modified Folch method (0.243pmol), indicating that modified SAPs method is more sensitive than the modified Folch method in positive-ion mode. For the possibility tests of the lipid extraction using the spin column with SAPs, lipid

extracts in another biological sample such as CSF were analyzed for lipid profiling. These results similarly corresponded with modified Folch methods.

Conclusions

Generally, conventional lipid extraction demands multiple steps, including drying, incubation, partitioning, and centrifugation, which results in tedious labor and low performance. However, the extraction of lipids using the spin column-based SAP is much simpler and faster because it eliminates the complicated processes. Therefore, we firmly conclude that the modified SAPs method is a very easy, fast, and quantitative method that is especially important due to its use of small amounts of samples from complex biological specimens, such as plasma/serum for clinical use.

Novel Aspect

A newly high-speed and quantitative lipid extraction method was developed; this method has extensive application prosperity in the clinical field.

References

1. J. F. Brouwers, E. A. Vernooij, A. G. Tielens, and L. M. van Golde, *J. Lipid Res.*, 40, 164 (1999)
2. Bang, G., Kim, Y. H., Yoon, J., Yu, Y. J., Chung, S., and Kim, J. A, *Anal. Chem.*, 89, 24 (2017)

ThP-135 / LIPIDOMIC APPROACH FOR ALGAL OIL CHARACTERIZATION BY FTICR MASS SPECTROMETRY

Stefano Chiaberge - Roberta Miglio - Davide Deriu - Sara Perucchini - Andrea Siviero

eni s.p.a., Renewable Energy & Environmental R&D, Novara, Italy

Keywords: Microalgae, Bio-Fuels, FTICR-MS, Lipidomics, Lipid Extraction

Introduction

Microalgae *Nannochloropsis salina* provides high lipid yields and, therefore, it can be used as feedstock for bio-fuel production. The solvent extraction of apolar lipids (mainly triglycerides-TAGs) from microalgae must be optimized for the transformation of these lipids into a bio-fuel through hydrotreating. The entire process optimization needs a detailed characterization of the microalgae extract(algal oil). A lipidomic[1] approach must be developed.

Methods

Algal oils were separated by solid phase extraction (SPE) using an amino based cartridge, obtaining two fractions: neutral compounds were eluted using a solvent mixture of chloroform and 2-propanol, polar compounds were obtained eluting with methanol. The fraction were then analyzed by flow injection with Electrospray ion source (ESI+) on a 7T-FTICR Mass spectrometer. The mass spectra were then elaborated[2]and the main peaks were attributed to a molecular formula.

Results

The algal oil was separated by SPE in two fractions which were dried and weighted. Starting from about 12mg of algae ethanol extract, about 9.1mg were collected as neutral lipids (74%), 1.6mg as polar lipids (13%). Therefore, a fraction of about 13% was lost in the sample loading step or bonded to the solid phase of the cartridge. The Apolar fraction was first analyzed by GCHT-FID for the quantitation of the main apolar lipid compounds. The starting extract and the two fractions were then analyzed by direct flow ESI+ FTICR MS for the semiquantitative determination of the main species (TAGs, Phospholipids and Glycolipids) mainly ionized as adducts with sodium. Thanks to the high resolution and mass accuracy of the FTICR mass spectrometer, the main species can be separated and most of the peaks can be attributed to a molecular formula. Tandem mass spectrometry allowed fatty acid identification for the main species.

The lipidomic approach developed provides a molecular fingerprint of the algal oil that is extremely useful for the optimization of the whole process.

Conclusions

A detailed characterization of algal oil was performed employing a lipidomic approach. Even if the SPE fractionation did not provide a total recovery, it allowed to obtain a preliminary quantification of apolar lipids (mainly TAGs) in the algal oil, by simple weighting of the extract. Finally, the detailed molecular characterization of lipid classes allowed the determination of the more suitable process conditions for the valorization of TAGs into a fuel-like mixture.

Novel Aspect

A novel analytical approach was developed for the direct quantification of apolar and polar lipids obtained by microalgae.

References

1. Liu P., Corilo Y. E., Marshall A. G., Energy & Fuels 30, 8270-8276 (2016)
2. Chiaberge S. et al. ChemSusChem 5, 1-9 (2012)

ThP-136 / DUAL MASS SPECTROMETRY AS A TOOL TO IMPROVE ANNOTATION AND QUANTIFICATION IN TARGETED PLASMA LIPIDOMICS

Liang Gao ⁽¹⁾ - **Amaury Cazenave-Gassiot** ⁽²⁾ - **Bo Burla** ⁽¹⁾ - **Markus Wenk** ⁽²⁾ - **Federico Torta** ⁽²⁾

National University of Singapore, Life Sciences Institute, Singapore, Singapore (1) - National University of Singapore, Department of Biochemistry, Singapore, Singapore (2)

Keywords: Dual MS, lipidomics, glycerophosphocholines, sphingomyelins

Introduction

Lipidomics is rapidly developing as one of the most important areas in biomedical research [1-2]. New technologies and analytical approaches are advancing a field that is defining new guidelines and standards. High quality data, based on reliable quantification and clear identification of the reported lipid species, are required for the clinical translation of human plasma lipidomic studies [3-4].

Methods

As a tool to improve the quality of targeted lipidomics studies, we established a Dual MS platform by simultaneously coupling a reversed-phase liquid chromatography (RPLC) separation to a triple quadrupole (QqQ) and a quadrupole-time of flight (Q-ToF) mass spectrometers. In one single experiment, this platform allows to correctly identify, by high resolution MS/MS, the peaks that are quantified by MRM.

Results

Glycerophosphocholines (GPCs) and sphingomyelins (SMs) play important biological roles in various cellular functions and physiological processes and are highly abundant in plasma. Quantification of GPCs and SMs can be efficiently performed on QqQ mass spectrometers in targeted multiple reaction monitoring (MRM) mode. A series of issues encountered when aiming at the identification and accurate quantification of GPC and SM in MRM were solved with the Dual MS platform, including (i) resolving peaks of GPC polyunsaturated species, (ii) discriminating between plasmanyl-, plasmenyl- and odd chain GPCs and (iii) resolving the isotopic overlap between GPC and SM co-eluting species.

Conclusions

Our platform provides a higher level of confidence in the quantification and identification of

phosphocholine-containing lipids in human plasma. Our work is aligned with several new initiatives in the lipidomic field with the aim to establish better methodologies and standard procedures to generate high quality lipidomic data.

Novel Aspect

This new methodology might have a great potential in comprehensive studies that want to achieve an acceptable specificity in the quantification of lipids.

References

1. Han X., Nat. Rev. Endocrinol. 12, 668-679 (2016).
2. McFadyen J. D., Meikle P. J., Peter K., Eur. Heart J. 38, 2006-2008 (2017).
3. Holčapek M., Liebisch G., Ekroos K., Anal. Chem. 90, 4249-4257 (2018).
4. Liebisch G., Ekroos K., Hermansson M., Ejsing C. S., Biochim. Biophys. Acta 1862, 747-751 (2017).

ThP-137 / NEW PEAK ALIGNMENT SOFTWARE: METABOALIGN IMPROVING ISSUES OF LC-MS/MS LIPIDOMICS ENABLES HIGHTHROUGHPUT DATA ANALYSIS FOR THOUSANDS OF LIPIDS IN BIOLOGICAL MATERIALS.

Akio Hayashi⁽¹⁾ - Takashi Matsuzaka⁽²⁾ - Hitoshi Shimano⁽²⁾

ONO Pharmaceutical Co.,LTD., Minase Research Institute, Mishima-gun, Japan (1) - University of Tsukuba, Faculty of Medicine, Tsukuba, Japan (2)

Keywords: ELOVL6, insulin resistance, targeted lipidomics, software, automatic-annotation

Introduction

Up to now many quantification softwares for targeted LC-MS/MS analysis using triple quadrupole mass spectrometer have been developed, but any of them could not be completely sufficient in peak alignment because of the difficulty to correct for any drift in retention time. Furthermore, it is difficult to identify the fatty acid composition of phospholipids (PLs) by simultaneous MS/MS assay using a thousand of polar head group MRM transitions because of the limitation of MS/MS cycle time in one assay. Therefore, we developed a new software for targeted lipidomics (TL) to improve these issues. Then we applied this software to lipidomics in liver of liver specific ELOVL6 knockout (LKO) mice, and found a lipid which was correlated with insulin sensitivity scores.

Methods

The liver tissues of LKO mice were homogenized in 1% acetic acid methanol solution on ice cold water. Then the whole lipid extraction and TL were performed by the method previously reported. TL and non-TL raw data were obtained by 4000QTRAP (Sciex) and QExactive HF (Thermo Fisher Scientific). The lipid identification data set including the information of PLs fatty acid analysis were made by Lipid Search (LS) software (Mitsui Knowledge Industry Co., Ltd.) from non-TL data. TL raw data analysis was performed by newly developed lipidomics software: MetaboAlign beta version with automatic peak annotation by means of merging with the LS identification data set.

Results

MetaboAlign was allow to align peaks independent of peak shape and peak width. All Lyso-PL regio isomers of sn-1 form and sn-2 form could be separated, as for example LPCacyl (18:1) of sn-1 and sn-2 were separated accurately and partially LPCacyl (cis-Vaccenyl) and LPCacyl (Oleyl). LPCacyl (cis-Vaccenyl, sn-2) was correlated well with LPCacyl (16:1, sn-2). *N*-Acylethanolamine (18:1) (AEA) was separated into two peaks which were considered AEA (Oleyl) and AEA (cis-Vaccenyl). Liver PCs with a variety of fatty acids were separated by high resolution C18 reverse phase column, aligned, quantified, and identified about fifty PC

species automatically by MetaboAlign. ELOVL6 is an enzyme which can elongate fatty acid chain from C16 to C18. LKO mice showed lower levels of C18 fatty-acyl PCs in liver than WT mice under high sucrose diet condition. Then we found a lipid_Y which was correlated with insulin sensitivity scores which might be involved in improvement of insulin resistance.

Conclusions

We successfully developed a new software MetaboAlign beta version for LC-MS/MS lipidomics to enable to align, quantify and identify thousands of lipid peaks. Both four isomers of LPCacyl (18:1) and two isomers of AEA(18:1) were aligned and quantified adequately. We found a lipid_Y which was correlated well with insulin sensitivity.

Novel Aspect

Successful development of a new software MetaboAlign beta version for LC-MS/MS lipidomics to enable to align, quantify and identify thousands of lipid peaks.

Separation of four isomers of LPCacyl (18:1) and two isomers of AEA (18:1).

Finding of a lipid_Y which was correlated with insulin sensitivity scores which might be involved in improvement of insulin resistance.

References

1. H. Saga, A. Hayashi, et al., PLoS ONE, 9, e93230(2014).
2. H. Muranaka, A. Hayashi, et al., Oncogenesis, 6, e350 (2017).

For information please contact: a.hayashi@ono.co.jp

ThP-138 / STRUCTURAL ANALYSIS OF LIPIDS FROM CELLS OF BACTERIA OF THE GOOD, THE BAD AND OF THE UGLY

Fong-Fu Hsu⁽¹⁾ - **Cheryl Frankfater**⁽¹⁾ - **William Henson**⁽²⁾ - **Albert Haas**⁽³⁾ - **John Turk**⁽¹⁾ - **Georgiana Purdy**⁽⁴⁾ - **R Tatituri**⁽⁵⁾

Washington University School of Medicine, Department of Medicine, St. Louis (1) - Washington University, Department of Energy, Environmental & Chemical Engineering, St. Louis (2) - University of Bonn, Molecular Microbiology, Bonn (3) - Oregon Health & Science University, Molecular Microbiology & Immunology, Portland (4) - Harvard Medical School, Division of Rheumatology, Immunology, and Allergy, Boston (5)

Keywords: lipids, bacteria, linear ion trap, high resolution mass spectrometry.

Introduction

Bacteria, regardless of whether they are beneficial or pathogenic all contain a large diversity of lipids with various structures and functions in the cells; and the structures are complex and difficult to define. Herein, we report LIT MSn method to identify the lipid structures from *R. opacus*, *Cordyceps sinensis*, *P. gingivalis*, *L. monocytogenes* and *Corynebacterium stratum*, which contain unique lipids different from those in eukaryotic and other prokaryotic cells.

Methods

Cells were homogenized and lipids were extracted using a modified Bligh & Dyer method. Preparative HPLC experiments were carried out using a Thermo TSQ Vantage MS with Accela UPLC. Separation of lipids was achieved by a Supelco C-8 column. High resolution LIT MSn experiments were performed on a LTQ Orbitrap Velos. Peptidolipids collected from HPLC separation were dried, re-dissolved in CHCl₃/CH₃OH, and continually infused into the ESI source.

Results

In addition to the lipids commonly seen for other bacteria, *R. Opacus* contains two novel families of

peptidolipids; *Cordyceps sinensis* contains extremely rich linoleic acid containing phospholipids; *P. Gingivalis* and *L. monocytogenes* contain odd chain methyl branch iso-fatty acid containing phospholipids, with the former also contains dihydroceramides, including EPC, IPC, GPC, and GS-peptidolipid, and the latter also contains plasmalogen PG and lysyl-CL; *C. stratum* contains a novel glucuronosyl diacylglycerol family. All the fatty acyl substituents in EPC, IPC, and GPC belong to the β -OH fatty acid family, and the structure is readily identifiable by CID MS2.

Conclusions

We found three new lipid families (peptidolipid, glycine-lipid, 18:1/16:0-glucuronosyl diacylglycerol) that have not been previously reported for these bacteria. ESI LIT CID/HCD MSn with high resolution mass spectrometry revealed the complexity of microbial lipid structures, including many isomers. The MSn capability of LIT adds additional dimensions in segregating ions, thereby, accelerating and simplifying data interpretation and structural characterization.

Novel Aspect

LIT MSn with high resolution MS unravels the structural details of novel lipids including many minor isomers that would be very difficult to define using other analytical methods.

ThP-139 / JOINT APPLICATION OF ESI AND MALDI LIPIDOMIC PROFILES FOR CLASSIFICATION OF BRAIN TUMORS

Daniil Ivanov⁽¹⁾ - *Anatoly Sorokin*⁽¹⁾ - *Evgeny Zhvansky*⁽¹⁾ - *Vsevolod Shurkhay*⁽²⁾ - *Stanislav Pekov*⁽¹⁾ - *Igor Popov*⁽¹⁾ - *Evgeny Nikolaev*⁽³⁾

Moscow Institute of Physics and Technology (State University), Moscow, Russian Federation (1) - Federal State Autonomous Institution "N. N. Burdenko National Scientific and Practical Center for Neurosurgery" of the Ministry of Healthcare of the Ru, Moscow, Russian Federation (2) - Skolkovo Institute of Science and Technology, Skolkovo, Russian Federation (3)

Keywords: brain tumors, lipidomics, molecular profiling, classification

Introduction

Analysis of molecular profiles is one of the promising approaches to determine malignancy of brain tissues. Different ionisation methods could be used for molecular profiling: online extraction with ESI is better for fast intraoperative analysis, while MALDI more suitable for histological investigation. The purpose of this study is to quantify ability to construct a joint classifier of brain tumors, working with MALDI and direct spray MS profiles simultaneously.

Methods

Set of molecular profiles of different brain tumor tissues was obtained by online extraction with ESI and MALDI Imaging experiments. Discriminating features for the malignancy level of tumor determination were identified individually for each ionisation method [1], and two brain tissue classifiers were created, followed by the comparison of selected features for each classifier.

Results

The molecular profiles that were obtained by direct ionisation and MALDI show the high level of correlation between each other for the same type of tissue especially in a mass range of common lipids. MALDI brain tissue classifier and direct spray classifiers were able to determine features to quantify the ratio of malignancy, and a number of these features were presented in both classifiers. Since we use high-resolution mass analysers for measuring direct-spray molecular profiles, for some of the common features we predict the molecular origin of these ions.

Conclusions

It was demonstrated that the classifiers build on the molecular profiles from direct-extraction MS and MALDI MS contain signals of the same compounds so both types of data could be used for the creation of joint classifier. The analysis of features, used in the classifiers, will allow optimising the parameters of experiments to improve the efficiency of classification and visualization of the spatial structure of the samples.

Novel Aspect

This work shows principal ability to use heterogeneous data to construct the classifier which works with online extraction and MALDI imaging experiments together.

References

1. Sorokin A., Zhvansky E., Shurkhay V., Bocharov K., Popov I., Levin N., Zubtsov D., Bormotov D., Kostyukevich Y., Potapov A., Nikolaev E., *European Journal of Mass Spectrometry*, 23, 237-241 (2017).

ThP-140 / DEVELOPMENT AND APPLICATION OF AN ENHANCED LIPIDOMIC PROFILING METHOD USING LC-MS

Malak Jaber - David Barrett - Dong-Hyun Kim

University of Nottingham, University of Nottingham, Advanced Materials and Healthcare Technologies, Nottingham, United Kingdom

Keywords: Lipidomics, Variations, ^{13}C in vivo labeling, Normalisation

Introduction

Quantification is a critical task in comprehensive lipidomics studies. Liquid chromatography-mass spectrometry (LC-MS)-based technique can provide a valuable tool for monitoring lipid changes as well as to identify specific biomarkers. A key limitation of this technique is the nonlinear responses due to stability, matrix effect, instrumental variations and others. However, those can be overcome by our novel method using ^{13}C -labelled organism extracts (1).

Methods

The lipidome profiles of five different microorganisms were compared to obtain the most coverage of human lipids leading to comprehensive compound-specific normalisation in quantitative mass spectrometry assays. An in vivo labelling strategy was used to produce a wide range of isotopically labelled lipid species from the selected microorganism grown on uniformly ^{13}C -labelled glucose, which have similar physico-chemical properties of unlabelled lipids (2, 3).

Results

Comparison of the lipid profiles of five different species including *E. coli*, *Spirulina*, *S. cerevisiae* CEN. PK 113-7D, *S. cerevisiae* BY4741, and *P. pastoris* NCYC 175 showed that there were marked differences in the relative amount of the identified lipid species as well to their labelling pattern. *P. pastoris* provides the highest number of common lipids with human metabolome database with good labelling efficacy. Therefore, *P. pastoris* was selected as a source of isotopically labelled internal standards (IS).

Then, the lipidomic extract of *P. pastoris* used to normalise variations introduced during samples preparation, storage and LC-MS analysis. The labelled yeast extract was successfully used to correct variations in LC-MS response in long analysis time (230 injections over 3.5 days) and reduction in CV% of the identified lipids was achieved (from 8.4 to 4.9). Also, the extract was able to reduce batch to batch variations successfully (from 17.2 to 6.2) without the need for quality control samples (QC) and thus decrease the analysis time.

Conclusions

P. pastoris proved to be the ideal source for production of labelled IS with more coverage for higher eukaryotic lipids. Normalisation using the ^{13}C -labelled extract could successfully correct for instrumental variation and batch to batch variations without the need for QC samples. These results show the high potential of the use of in vivo labelling strategy in quantitative measurement of a wide range of lipids especially in more complex matrices.

Novel Aspect

To use labelled organism extracts as a source of IS to develop more reproducible, quantitative untargeted lipidomic analysis that can be applied to a variety of clinical studies

References

1. Lei Z, Huhman DV, Sumner LW, Journal of Biological Chemistry, 286, 25435-25442 (2011).
2. Rampler E, Criscuolo A, Zeller M, El Abiead Y, Schoeny H, Hermann G, et al., Analytical Chemistry, 90, 6494-6501 (2018).
3. Rampler E, Coman C, Hermann G, Sickmann A, Ahrends R, Koellensperger G., Analyst, 142, 1891-1899 (2017).

ThP-141 / UPLC-QQQ/MS-BASED LIPIDOMICS ANALYSIS OF HIGH GLUCOSE INDUCED HEPG2 CELLS TREATED WITH HERBAL MEDICINE EXTRACTS

Jae Won Lee - Dojin Choi - Bo Ram Choi - Youngseob Lee - Geumsoog Kim - Daeyoung Lee

National Institute of Horticultural and Herbal Science, Department of Herbal Crop Research, Eumseong, Republic of Korea

Keywords: Lipidomics, HepG2 cells, oxidative stress, LC/MS, multivariate analysis

Introduction

Liver is the primary organ for glucose metabolism and regulation; therefore it is interesting to elucidate the effects of hyperglycemia in cultured liver cells. HepG2 cell, the human hepatoma cell line, has been used extensively to study hyperglycemia in vitro as evidenced by several reports [1,2]. HepG2 cells treated with glucose exhibited an upregulation of the transcription of human apolipoprotein A-II gene which causes a rise in plasma triglyceride level and glucose intolerance, resulting in hyperglycemia [3]. High glucose exposed HepG2 cells exhibited inhibition of 50AMP-activated protein kinase (AMPK) which plays a role in cellular energy homeostasis, and increased accumulation of lipids which were prevented either by metformin or by expressing the constitutively active AMPKa [4]. The high concentration of glucose also leads to oxidative stress, which plays an important role in the glucose-induced cellular dysfunction in diabetes complications [5]. Oxidative stress is defined as an imbalance between free radical production and antioxidant defense in a biologic system. It has been reported that irregular cellular metabolism in diabetes leads to generation of free oxygen radicals and reduction of antioxidant capacity [6].

In cells, lipids have functions as structural membrane components, sources of energy, and signaling molecules. As hormones and cellular messengers, lipids are particularly involved in many intracellular and intercellular signaling processes [7]. Furthermore, lipids adjust membrane fluidity and support diverse cellular processes [8]. Thus, it is critical to determine the cellular lipid metabolism correlated with oxidative stress.

Recently, liquid chromatography coupled to tandem mass spectrometry (LC-MS/MS) has been used as an effective tool to analyze various lipids from biological samples [9]. An ultraperformance LC (UPLC) enables the rapid and effective separation of individual lipid species [10]. Multiple reaction monitoring (MRM) based on triple quadrupole (QqQ)/MS is also applicable for the sensitive analysis of target compounds [11].

Optimal lipidomics based on UPLC-QqQ/MS has been widely applied to shed light on many biological events [12]. In this study, UPLC-QqQ/MS-based lipidomics was used to elucidate the altered lipid metabolism of HepG2 cells affected by the oxidative stress. We also treated anti-oxidative compound and extracts to the high glucose induced HepG2 cells and assessed the altered levels of lipids. The detailed characterization of lipids in HepG2 cells having oxidative stresses may provide a good understanding of the lipid metabolism correlated with oxidative stress.

Methods

MTT and Lactate dehydrogenase (LDH) assay

Cell viability was determined by the quantitative colorimetric assay using MTT reagent. The cytotoxicity was also measured by lactate dehydrogenase (LDH) assay using commercial LDH detection kit for detection the leakage of soluble cytoplasmic LDH into the extracellular medium.

Lipid extraction

First, the cell pellet was added to 990 μL chloroform/methanol (1:2, v/v) and 10 μL 1 $\mu\text{g}/\text{mL}$ of lipid standards. The sample was incubated for 10 min on ice with vortexing for 30 s. After centrifugation (13,800 \times g, 2 min at 4°C), 950 μL supernatant were transferred to a new Eppendorf tube. Second, the remaining pellet was resuspended in 750 μL chloroform/methanol/37% (1N) HCl (40:80:1, v/v/v) and incubated for 15 min at room temperature with vortexing. The tube was transferred to ice, and 250 μL cold chloroform and 450 μL cold 0.1 M HCl were added followed by 1 min of vortexing and centrifugation (6,500 \times g, 2 min at 4°C). The lower organic phase was added to the tube, and subjected into LC/MS.

UPLC-QqQ/MS analysis

Lipid profiling was performed using the 6490 Accurate-Mass Triple Quadrupole LC-MS coupled to a 1200 series HPLC system (Agilent Technologies, Wilmington, DE, USA) with a Hypersil GOLD column (2.1 \times 100 mm ID; 1.9 μm , Thermo Science). Solvent A consisted of an acetonitrile–methanol–water mixture (19:19:2) with 20 mmol/L ammonium formate and 0.1% (v/v) formic acid, and solvent B consisted of 2-propanol with 20 mmol/L ammonium formate and 0.1% (v/v) formic acid. The flow rate was 250 $\mu\text{L}/\text{min}$, and total run time was 30 min for each analysis. MRM conditions including transition and MS/MS collision energy were optimized to analyze target lipids in individual samples.

Data processing and multivariate analysis

LC/MS data were obtained and processed using Agilent Mass Hunter Workstation Data Acquisition software. The MRM data of target lipids, including m/z of precursor ions, m/z of product ions, and retention time (RT), were exported using Qualitative Analysis B.06.00 software (Agilent Technologies, USA). Next, the Skyline software package (MacCoss Laboratory, University of Washington, USA) was used as an in-house database to determine the peak area of each assigned lipid from replicate raw data. The extracted areas of lipid peaks were normalized to the appropriate IS. Multivariate analyses such as principal component analysis (PCA), the sparse projection to latent structure discriminant analysis (sPLS-DA), and analysis of variance (ANOVA) were performed using the MetaboAnalyst website.

Results

Effects of ALM16 on high glucose-induced cell death and cytotoxicity in HepG2 cells

To determine the protective effects of the ALM16 on high glucose-induced cell death and cytotoxicity, HepG2 cells were pre-treated with the ALM16 for 2 h, and then treated with high glucose (30 mM) for 24 h. The ALM16 (10-200 $\mu\text{g}/\text{mL}$) did not show significant cytotoxic effects. While treatment of the HepG2 cells with high glucose (30 mM) significantly reduced the cell viability to approximately 64.0% compared with the untreated control cells (100%), pre-treatment with the ALM16 (50, 100, and 200 $\mu\text{g}/\text{mL}$) for 2 hours before exposure to high glucose (30 mM) dose-dependently reduced the high glucose-induced cell death to 70.2%, 73.3%, and 78.5%, respectively. In addition, pre-treatment with silymarin (50 $\mu\text{g}/\text{mL}$), as a positive control, reduced the cell death to approximately 83.0%. The effect of the ALM16 on the high glucose-induced hepatotoxicity in HepG2 cells was confirmed by the LDH assay. Whereas an approximately 2.89-fold higher LDH activity was observed from the high glucose-treated cells compared with the untreated control cells (1.00 fold), HepG2 cells pre-treated with the various concentrations of ALM16 were reduced in a dose-dependently manner by 2.78-, 2.40-, and 2.23-fold increases in LDH activity, respectively. The cells

pre-treated with silymarin (50 $\mu\text{g/ml}$) were also reduced LDH activity to 2.53-fold increase.

Lipid profiling of high glucose induced HepG2 cells treated with the mixture of *Astragalus membranaceus* and *Lithospermum erythrorhizon* extracts

Next, six samples of HepG2 cells (normal control (Con), Con + glucose 30 $\mu\text{g/ml}$ (Glu30), Glu30 + silymarin 50 $\mu\text{g/ml}$ (Sil50), Glu30 + ALM16 50 $\mu\text{g/ml}$ (Alm50), Glu30 + ALM16 100 $\mu\text{g/ml}$ (Alm100), Glu30 + ALM16 200 $\mu\text{g/ml}$ (Alm200)) were subjected into lipidomic analysis. For the profiling of various lipids from these HepG2 cells, we used a previously constructed method using UPLC-QqQ/MS [13]. Lipid standards were used to optimize the MRM conditions for the analysis of TG, DG, MG, ChE, cholesterol, PC, PE, PS, PG, PI, PA, LPC, LPE, LPS, LPG, LPI, LPA, SM, Cer, dCer, So, Sa, Cer1P, So1P, and Sa1P. In particular, TMSD was used to methylate several highly acidic lipids, such as PS, PI, PA, LPS, LPI, LPA, Cer1P, So1P, and Sa1P. By using an octadecylsilyl silica column, various lipid species were separated according to their total carbon chain length (Cn) and total degree of unsaturation (Un). In the MRM data, the RT of each lipid was increased with a higher Cn and a lower Un. In this study, the skyline software was used to assign each lipid peak in the basis of RT and to calculate its peak area in the MRM data. Finally, in the HepG2 cell, a total of 336 lipid species were identified as follows: 23 TGs, 37 DGs, 8 MGs, 16 ChEs, cholesterol, 62 PCs, 28 PEs, 15 PSs, 5 PGs, 16 PIs, 34 PAs, 17 LPCs, 13 LPEs, 5 LPSs, 7 LPGs, 4 LPIs, 15 LPAs, 11 SMs, 11 Cers, 2 dCers, So, Sa, 2 Cer1Ps, So1P, and Sa1P. Furthermore, in the data processing of each lipid class, the calculated peak area of individual lipid species was normalized by the peak area of the IS. And then, the normalized datasets were subjected to statistical analysis.

Statistical analysis of lipid datasets obtained by the profiling of six HepG2 cells

The six datasets of lipids analyzed from Con, Glu30, Sil50, Alm50, Alm100, and Alm200 were subjected into multivariate analyses. Assays were carried out in triplicate ($n = 3$) for each sample. First, we applied the multiple pattern recognition methods such as PCA and sPLS-DA for visualizing the general clustering trends among the six groups. In the PCA score plot, the six groups were not sufficiently distinguished. For well discriminating the six groups, we also applied the sPLS-DA that reduce the number of variables in the data to produce robust and easy-to-interpret model. In the score plot, points of Con were firstly separated with those of Glu30. Furthermore, other four groups including Sil50, Alm50, Alm100, and Alm200 were scattered between Con and Glu30 in the basis of Component 1. These indicated that cell lipid profiles differ depending on normal control, high glucose induced cells, and high glucose induced cells treated with silymarin or ALM16. In the sPLS-DA score plot, Sil50s were getting closer to Cons. This indicated that the lipid profiles of Sil50 are similar with those of Con comparing to Glu30. Furthermore, the points of Alm200, Alm100, and Alm50 were in order getting closer to those of Con. When the treatment amounts of ALM16 were increased, their points were closer to Cons than others. This indicated that ALM16 treatment is effective to suppress the cellular oxidation and recover the lipid metabolism. In the loading, we also presented the variable importance in projection (VIP) plot of 25 lipids (VIP scores Top 25) that were differently regulated among the six samples of HepG2 cells. Next, ANOVA test was performed to find the differentially regulated lipids (DRLs), which were statistically different lipids among the six groups. As a result, the p-values of 142 lipids were below 0.05 by using the Student's t-test, and we selected several lipids that have the lowest p-values to find the alteration of their relative levels in the six groups.

Conclusions

In this study, UPLC-QqQ/MS-based lipidomics was used to elucidate the altered lipid metabolism of HepG2 cells affected by the oxidative stress. In the cell experiments, we confirmed that high glucose treatment causes the cellular oxidation. Then, ROS may give effects on the metabolism of various metabolites including lipids. Silymarin is a well-known drug for the anti-oxidation activity. Thus, the treatment of silymarin may suppress the oxidation of high glucose induced cells. Alm treatment also showed the biological activity to suppress the oxidative stress. Thus, the treatment of silymarin or Alm may suppress the severe alteration of lipid metabolism in high glucose induced HepG2 cells. To detail characterize the lipid alterations of HepG2 cells affected by oxidative stress and treatment of silymarin or Alm, six samples of Con, Glu30, Sil50, Alm50, Alm100, and Alm200 were subjected into lipidomic analysis. Furthermore, multivariate analyses such as PCA, sPLS-DA, and ANOVA were applied to characterize the lipid alterations in these six groups. As a result, sPLS-DA well differentiated the lipid profiles of these six groups. In the ANOVA

tests, 142 DRLs were identified to find the alteration of their relative levels in the six groups. Finally, we assessed the altered levels of lipids by oxidative stress, and suggested these lipids as the biomarker candidates associated with oxidative stress. Furthermore, we assessed the lipid alteration in the high glucose induced HepG2 cells treated by silymarin or Alm. When silymarin or Alm suppresses the damage of oxidative stress, it can decrease the degree of lipid alteration occurred by oxidative stress. With these results, it was possible to more understand how lipid metabolism of HepG2 cells was affected by oxidative stress. In future aspect, identified DRLs can be used as biomarker candidates to evaluate the oxidation of liver cells.

Novel Aspect

This study firstly attempted to perform the LC/MS-based lipidomics analysis of high glucose induced HepG2 cells which were damaged by oxidative stress. Various lipids including glycerol lipids (*i.e.*, TG, DG, MG), sterol lipids (*i.e.*, ChE, cholesterol), phospholipids (*i.e.*, PC, PE, PS, PG, PI, PA), lysophospholipids (*i.e.*, LPC, LPE, LPS, LPG, LPI, LPA), and sphingolipids (*i.e.*, SM, Cer, dCer, So, Sa, Cer1P, So1P, Sa1P) were profiled in the HepG2 cells, and their levels were assessed in six samples of Con, Glu30, Sil50, Alm50, Alm100, and Alm200. We also found how the pretreatment with silymarin or Alm before exposure to high glucose dose dependently attenuated the high glucose induced hepatotoxicity and altered lipid metabolism. The DRLs identified by multivariate analyses can be used as biomarker candidates to evaluate the oxidation of liver cells.

References

1. Dey A., *Emergency Medicine: Open Access*, 3, 5 (2013)
2. Donato M. T., Tolosa L., Gómez-Lechón M. J., *methods in Molecular Biology*, 1250, 77-93 (2015).
3. Sauvaget D., Chauffeton V., Dugué-Pujol S., Kalopissis A. D., Guillet-Deniau I., Foufelle F., Chambaz J., Leturque A., Cardot P., Ribeiro A., *Diabetes*, 53, 672-8 (2004).
4. Zang M., Zuccollo A., Hou X., Nagata D., Walsh K., Herscovitz H., Brecher P., Ruderman N. B., Cohen R. A., *Journal of Biological Chemistry*, 279, 47898-905 (2004).
5. Li Q., Lin Y., Wang S., Zhang L., Guo L., *Scientific Report* 7, 8008 (2017).
6. Matough F. A., Budin S. B., Hamid, Z. A., Alwahaibi N., Mohamed J., *Sultan Qaboos University Medical Journal*, 12, 5-18 (2012).
7. Higdon A., Diers A. R., Oh J. Y., Landar A., Darley-Usmar V. M., *Biochemical Journal*, 442, 453-464 (2012).
8. Csoboz B., Balogh G. E., Kusz E., Gombos I., Peter M., Crul T., Gungor B., Haracska L., Bogdanovics G., Torok Z., *International Journal of Hyperthermia*, 29, 491-499 (2013).
9. Scherer M., Leuthäuser-Jaschinski K., Ecker J., Schmitz G., Liebisch G., *Journal of Lipid Research*, 51, 2001-2011 (2010).
10. Gao X., Zhang Q., Meng D., Isaac G., Zhao R., Fillmore T. L., Chu R. K., Zhou J., Tang K., Hu Z., *Analytical and Bioanalytical Chemistry*, 402, 2923-2933 (2012).
11. Kitteringham N. R., Jenkins R. E., Lane C. S., Elliott V. L., Park B. K., *Journal of Chromatography B*, 877, 1229-1239 (2009).
12. Gaudin M., Panchal M., Auzeil N., Duyckaerts C., Brunelle A., Laprévotte O., Touboul D., *Bioanalysis*, 4, 2153-2159 (2012).
13. Lee J. W., Mok H. J., Lee D. Y., Park S. C., Kim G. S., Lee S. E., Lee Y. S., Kim K. P., Kim H. D., *Journal of Proteome research*, 16, 1460-1469 (2017).

ThP-142 / CONSTRUCTION OF NOVEL LIPIDOMICS PLATFORM COMBINATED OF TARGETED AND NON-TARGETED ANALYSIS

Yasuto Yokoi⁽¹⁾ - Hideya Kuwabara⁽¹⁾ - Tadahiro Hoshino⁽¹⁾ - Akio Hayashi⁽²⁾

Mitsui Knowledge Industry Co., Ltd., Bioscience Department, Tokyo, Japan (1) - ONO Pharmaceutical Co.,

Ltd., Minase Research Institute, Osaka, Japan (2)

Keywords: lipidomics; informatics; quantitative analysis; RT correction

Introduction

Targeted lipidomics (TL) using selected ion monitoring (SRM) enable high sensitive quantitative analysis for certain lipids in biological sample. However, overlapped chromatographic peaks from isomer and fluctuations in retention time of chromatograms cause miss alignment among samples. Here we developed novel lipidomics platform named MetaboAlign which can solve these issues by combining TL with Non-targeted lipidomics on data dependent acquisition (DDA).

Methods

MetaboAlign, which is consisted of peak detection module, alignment module, and viewer was written in Java. SRM data for rat plasma samples (n10) measured by 4000 Qtrap (Sciex) were converted to MzML by ProteoWizard. DDA data were measured by QExactive HF (Thermo Fisher Scientific), which were used for lipid identification by LipidSearch4.2™ (Mitsui Knowledge Industry). Both results were imported to MetaboAlign, quantitative profiles of lipids were calculated.

Results

MetaboAlign performs quantitative analysis including peak detection, alignment, identification and calculation of concentration automatically. RT correction algorithm enables large number of peaks in each sample to align with high accuracy. Combination with non-targeted lipidomics (NTL) using high resolution mass (HRMS) can also assign lipid to each peak in distinction from its isomer. We applied this software to rat liver sample (n10). As a result, 861 kinds of lipid species were identified by NTL, 496 of which were acquired as lipid profiles by TL.

Conclusions

Novel lipidomics platform combined TL with NTL was developed. Rapid and accurate quantitative analysis of large sample sets was successfully achieved using the software.

Novel Aspect

Comprehensive and effective platform for large scale quantitative analysis using combination of targeted and non-targeted lipidomics is presented.

ThPS-S07 - Developments and applications in ion mobility MS

Chairs: Perdita Barran, Leopoldo Ceraulo

ThP-143 / A HIGH-RESOLUTION ION MOBILITY SPECTROMETER COUPLED TO A TIME-OF-FLIGHT MASS SPECTROMETER

Maria Allers - Laila Timoumi - Ansgar T. Kirk - Florian Schlottmann - Stefan Zimmermann

Leibniz Universität Hannover, Institute of Electrical Engineering and Measurement Technology, Department of Sensors and Measurement Technology, Hannover

Keywords: High-resolution ion mobility spectrometer; Gated IMS-MS interface

Introduction

Over the last years, we have developed compact high-resolution drift tube ion mobility spectrometers (IMS) combined with different ion sources reaching high resolving power up to $R=250$ [1,2]. High-resolution

IMS allow fast separation of ion species by size, shape and charge, yielding structural information complementary to molecular mass. To keep such high IMS resolving power when coupling the IMS to a time-of-flight mass spectrometer (TOF-MS), here a Bruker micrOTOF II, several aspects have to be considered, such as the ion transfer and IMS-MS synchronization.

Methods

Usually, when coupling IMS to TOF-MS data acquisition of both instruments is synchronized resulting in two-dimensional data containing ion mobility and mass spectra. However, due to the short temporal peak widths of just 60 μ s in our high-resolution IMS with $R=100$, this technique is not applicable due to significant peak broadening during the ion transfer into the MS and an insufficient data acquisition rate of the MS. Thus, a simple gated interface has been designed to transfer selected peaks of the ion mobility spectrum into the TOF-MS. The interface is realized by replacing the Faraday plate of the IMS with a Faraday grid that is shielded by two additional aperture grids. By directly inserting the interface into the capillary shield of the MS, the IMS can be operated at ambient pressure and easily adapted to the MS inlet capillary.

Results

Coupling our high-resolution IMS to a TOF-MS as described above does not affect IMS resolving power. The novel interface allows both, recording full ion mobility spectra and transferring selected ion peaks or time domains of the ion mobility spectrum into the MS to analyze selected ion species formed in the IMS. Furthermore, due to high grid transparency, ion losses in the interface are minimized. The investigation of the positive reactant ions, ketones and alcohols show the occurrence of ion cluster dissociation and ion fragmentation processes inside the transfer stage of the Bruker micrOTOF II. As proof-of-concept study, the presented IMS-MS is used to investigate the product ions of benzene, toluene and m-xylene in dry, purified air using a 3H source. Besides the well-known product ions M^+ and $M\cdot NO^+$, for high concentrations a dimer ion is also observed for benzene and toluene, consisting of two molecules and three further hydrogen atoms.

Conclusions

A simple interface to couple high-resolution IMS to APCI-TOF-MS while keeping full IMS resolving power is presented. Due to its fast gating capability this interface can be used to investigate ion species of a certain peak in the ion mobility spectrum and to explore the ionization processes in IMS. In particular, coupling our IMS to a Bruker micrOTOF II allows investigating the ion chemistry of different ion sources.

Novel Aspect

Simple interface for coupling high-resolution IMS operated at atmospheric pressure to TOF-MS, minimizing ion losses, recording full ion mobility spectra and enabling a gated transmission of ions into the MS.

References

1. A. T. Kirk et al., *Analyst*, 138, 5200–5207 (2013)
2. A. T. Kirk et al., *Int. J. Ion Mobil. Spec.*, 18, 17–22 (2015)

Acknowledgement

This work was funded by the Deutsche Forschungsgemeinschaft (DFG, German Research Foundation) – ZI 1288/8-1

ThP-144 / APPLICATION OF PASEF MS/MS SCANS TO MONOCLONAL ANTIBODY PEPTIDE MAPPING

Anjali Alving⁽¹⁾ - Stuart Pengelley⁽²⁾ - Guillaume Tremintin⁽³⁾ - Detlev Sukau⁽²⁾

Bruker Daltonics, Biopharma, Billerica (1) - Bruker Daltonics, Biopharma, Bremen (2) - Bruker Daltonics, Biopharma, San Jose (3)

Introduction

Peptide mapping is an essential method for evaluating monoclonal antibodies (mAb) heterogeneity and their response to stress conditions. The demands placed on this assay have drastically evolved from initially checking the primary sequence to the localization of post translational modification and the determination of low level heterogeneities such as amino acid misincorporation. Because of this, there is a lot of interest in developing platform methods that are faster or able to provide more information.

Parallel serial fragmentation (PASEF) scans use ion mobility as a method to accumulate and select precursors for data dependent scans. This makes it possible to separate peptides that do not separate in the chromatography dimension, leading to improved sensitivity and new opportunities to improve peptide maps.

Methods

The NISTmab reference material (SRM8671) was incubated under high pH for 24 and 48h. The stressed samples and a control were alkylated and digested with trypsin before UHPLC separation using gradient lengths from 5 to 45 min. The data acquisition was conducted on a timsTOF Pro (Bruker) using the PASEF scan mode.

The resulting data was processed with Byonic/Byologic (Protein Metrics) and the ability to detect deamidated peptides was systematically evaluated for the various gradient lengths.

Preliminary Data

Rapid gradient: Sequence coverage upward of 95% was achieved for the unstressed NISTmab sample with an HPLC gradient length of 5 min and detection using the PASEF scan mode. Despite the short run time, multiple deamidation sites were detected.

Detailed characterization: The traditional peptide map method using an extended 120 min gradient combined with the PASEF scan mode allows the detection of low level heterogeneities down to 1-10 ppm of the NISTmab such as beta-microglobulin.

Now ongoing: We are now running a systematic evaluation of the performance of the method with gradient lengths of 5, 15, 30 and 45 min, while optimizing of the mobility and MS/MS parameters to take in account the increasing amount of coeluting precursors as the HPLC gradient length is reduced.

Novel Aspect

Impact of PASEF scans on the ability to detect deamidation in peptide mapping

ThP-145 / A CLIP ON ION MOBILITY DEVICE. ENABLING GAS PHASE IONS SEPARATION AND COLLISION CROSS SECTION MEASUREMENTS TO ANY MASS SPECTROMETERS

Bruno Bellina⁽¹⁾ - ***Jan Commandeur***⁽²⁾ - ***Perdita Barran***⁽¹⁾

University of Manchester, Manchester Institute of Biotechnology, Manchester (1) - MS Vision, Almere (2)

Keywords: Ion mobility, Instrumentation, Retrofitted

Introduction

Ion Mobility is a powerful gas phase separation technique that is used to distinguish ions as a function of their sizes/shapes; leading to molecule identification directly from the measurement of their mobilities. When coupled to mass spectrometry, IM-MS technique brings in a new dimension and numerous applications. We are developing a clip-on device that would be adaptable to any mass spectrometer to enable ion mobility separation and CCS measurements.

Methods

The developed prototypes utilise a resistive glass tube to generate an absolute and linear electric field, needed for direct collision cross sections (CCS) measurements. The entrance and the exit caps have been

specially designed to minimize hardware modifications on the host MS instrument. A single gas inlet is needed along with a total of only eight simple electrical connectors.

Results

The first generation will be fitted to a WATERS QToF-Ultima. The only necessary modification will consist into shortening the second ion guide, prior to the quadrupole region; retaining the original footprint of the instrument. An ion funnel has been designed and mounted on the exit cap allowing efficient refocusing of the ion beam post drift. A second generation will embed an ion mobility gating region (necessary for ion pulsing) along with the ion funnel region. Both have been developed to minimize the imprint on the host system and will be mounted and part of the entrance and exit caps. SIMION simulations have been successfully run on various types of model systems ranging from Leucine Enkephalin (low mass singly charged) to Cytochrome C (High mass, multiple charges). The maximum theoretical resolving power at different drift voltages at ambient temperature is evaluated with a maximum operating voltage of 20 V/cm.

Conclusions

A portable ion mobility device has been developed to enable gas phase separation and CCS measurement on any mass spectrometer (Q-ToF, FT and Ion Trap). The device is designed to minimize its imprint on the host, embedding trapping, drifting and refocusing regions in a single assembly. The simulations undertaken show very promising results for a large range of ions.

Novel Aspect

The novelty resides in the device flexibility. The length of the ion mobility drift region is easily adjustable and the device can be fitted to any MS instruments.

ThP-146 / EXPLORING THE CHEMICAL COVERAGE AND STRUCTURAL ASSIGNMENT OF NATURAL PRODUCTS BY ION MOBILITY SPECTROMETRY COMBINED WITH ULTRA-PERFORMANCE LIQUID CHROMATOGRAPHY/ MASS SPECTROMETRY

Fausto Carnevale Neto^(1,2) - *Joseph M. Egan*⁽²⁾ - *Leticia Veras Costa-Lotufo*⁽³⁾ - *Norberto Peoporine Lopes*⁽¹⁾ - *Roger G. Linington*⁽²⁾

Faculdade de Ciências Farmacêuticas da Universidade de São Paulo, Ribeirão Preto, São Paulo, Brazil (1) - Department of Chemistry, Simon Fraser University, Burnaby, BC, Canada (2) - Departamento de Farmacologia, Instituto de Ciências Biomédicas, Universidade de São Paulo, São Paulo, Brazil (3)

Keywords: Ion mobility-mass spectrometry, Natural products, CID fragmentation, Collisional cross section.

Introduction

The major goal in Natural Products (NP) research is to interrogate complex biological extracts for the purposes of metabolic exploration and natural products discovery. In order to achieve this objective, NP studies rely on accurate and selective acquisition of all possible chemical information, which includes maximization of the number of detected metabolites and their correct molecular assignment [1-3].

Methods

In this work, we investigated the capabilities of liquid chromatography coupled to ion mobility high-resolution mass spectrometry in the context of natural products. In other words, we addressed the use of comprehensive four-dimensional data (chromatographic retention time, drift time, mass-to-charge and ion abundance) as a way of detect and elucidate NP in complex biological extracts at nanoscale.

Results

In order to demonstrate the performance of ion mobility for larger chemical coverage, we investigated the

behavior of distinct NPs (or NP-inspired molecules) on UHPLC within hybrid Travelling Wave Ion Mobility (TWIM) QTOF system, considering sensitivity and selectivity, LC-deconvolution and conformation space (correlation between m/z and drift time) filtering. The procedure focused on two types of MS/MS data acquisition; data-dependent acquisition (DDA) or auto-MS/MS, and data-independent acquisition (DIA) or all-fragment-ion technique. We evaluated the influence of the “precursor ion isolation width” on chemical coverage, using both approaches, and the impact of additional ion mobility dimension on molecular detection and spectral-quality information. The observations were further explored in a partially characterized sample from a bank of marine bacterial extracts, spiked with the mixture of STDs at the two ranges of concentrations.

Conclusions

The use of ion mobility (IMS) as orthogonal information on LC-MS method can improve both metabolite detection and chemical characterization for natural products by separating ions in the gas-phase and adding coarse-grained structural information based on shape and collision cross sectional area (CCS).

Novel Aspect

The application of IMS-MS can improve both metabolite detection and chemical characterization of complex samples in the field of Natural Products.

References

1. Ortmayr K., Causon T. J., Hann S., Koellensperger G. *TrAC Trends Anal. Chem.*, 82, 358–366 (2016).
2. Kurita K. L., Glassey E., Linington R. G. *Proc. Natl. Acad. Sci. USA*, 112, 11999–12004 (2015).
3. Aksenov A. A., da Silva R., Knight R., Lopes N. P., Dorrestein P. C. *Nat. Rev. Chem.* 1, 54 (2017).

ThP-147 / INFLUENCE OF PROTOMERS ON DISSOCIATION OF SEQUENCE-DEFINED OLIGO(ALKOXYAMINE AMIDE)S

Christophe Chendo⁽¹⁾ - *Salomé Poyer*⁽²⁾ - *Didier Siri*⁽²⁾ - *Jean-François Lutz*⁽³⁾ - *Laurence Charles*⁽²⁾

Fédération des Sciences Chimiques (FR1739), Aix Marseille Université, Marseille (1) - Institut de Chimie Radicalaire (UMR7273), Aix Marseille Université, Marseille (2) - Institut Charles Sadron, CNRS, Chimie Moléculaire de Précision (UPR22), Université de Stasbourg, Strasbourg (3)

Keywords: Protomers; dissociation mechanism; ion mobility; sequence-defined polymers; MS/MS

Introduction

Poly(alcoxyamine amide)s are sequence-defined synthetic polymers in which binary information can be written thanks to a 0/1 coding system based on the mass of two amide synthons that alternate with a nitroxide moiety [1]. They can be easily sequenced by MS/MS since they mainly dissociate via homolysis of all fragile C–ON bonds between a coding unit and a nitroxide moiety. Surprisingly, each of these cleavages does not produce the two expected complementary products.

Methods

Oligo(alcoxyamine amide) trimeric models were ionized in positive ion mode electrospray. Tandem mass spectrometry and ion mobility spectrometry (IMS) experiments were performed with a Synapt G2 HDMS instrument. Geometry optimization was done after DFT calculations using M06-2X functional with the 6-31G(d) basis set implemented in the Gaussian-16 program. The MOBCAL software was used to convert Cartesian coordinates in collision cross section (CCS) values.

Results

Production of only one fragment upon homolysis of one X-T bond in a protonated α -X-T-X-T-X- ω trimer (where X stands for the coding amide moiety and T the TEMPO nitroxide) strongly suggests that, although

not directly involved in the mechanism of this charge-remote reaction, the adducted proton has an influence on its occurrence [2]. Protonation of the nitrogen atom in nitroxide was demonstrated to induce an increase of the dissociation energy of the C–ON bond [3]. Accordingly, homolysis of the first X-T bond would occur only when the second T nitroxide was protonated, and vice versa. This assumption implies that each precursor ion would exist as two protomers. Ion mobility experiments indeed showed two distinct conformations for all protonated trimers but those with a 1 coded unit in the second monomer, which were shown to have the same CCS regardless of the proton location. However, although a single C–ON bond homolysis was expected for each protomer, MS/MS spectra recorded after IMS separation still showed homolysis of the two C–ON linkages.

Conclusions

Positive ion mode electrospray ionization was shown to produce different protomers of sequence-defined oligo(alkoxyamine amide)s. By combining ion mobility data and molecular modeling, conformation of these protomers was demonstrated to be sequence-dependent. However, the actual influence of proton location on the occurrence of alkoxyamine bond homolysis is still unclear and is currently investigated with theoretical calculations.

Novel Aspect

Protomer-specific conformations and their influence on dissociation of small sequence-controlled polymers

References

1. Roy R. K., Meszynska A., Laure C., Charles L., Verchin C., Lutz J. -F., Nat. Commun. 6, 7237 (2015).
2. Charles L., Laure C., Lutz J. F., Roy R. K., Macromolecules 48, 4319-4328 (2015).
3. Mazarin M., Girod M., Viel S., Phan T. N. T., Marque S. R. A., Humbel S., Charles L., Macromolecules 42, 1849-1859 (2009).

ThP-148 / TRAPPED ION MOBILITY SPECTROMETRY AS NOVEL SEPARATION TECHNIQUE FOR THE INVESTIGATION ON ISOMER FORMATION DURING THE SIMULATION OF PHASE I METABOLISM BY MEANS OF ELECTROCHEMISTRY/MASS SPECTROMETRY

Jens Fangmeyer - Simon Gereon Scheeren - Robin Schmid - Uwe Karst

University of Muenster, Institute of Inorganic and Analytical Chemistry, Münster

Keywords: Electrochemistry, Metabolism, Trapped Ion Mobility Spectrometry, Liquid chromatography, Hyphenation

Introduction

Electrochemistry/mass spectrometry (EC/MS) has evolved as a valuable and purely instrumental approach in mimicking the phase I metabolism of xenobiotics. Due to the formation of isomers after the electrochemical treatment, the development of a separation method in order to distinguish between them is required. For this purpose trapped ion mobility spectrometry (TIMS) enables a novel and faster approach for the separation of oxidation products.

Methods

The generation of metabolites was performed in an electrochemical cell equipped with a boron-doped diamond working electrode applying a potential ramp. The EC cell was coupled online to ESI-HR-MS. The results were plotted in three-dimensional mass voltammograms. Mobile voltammograms for each m/z were recorded and visualized as a contour plot, which shows the signal intensity versus the applied potential and the inverted reduced mobility ($1/k_0$) of the analyte ions.

Results

With the recorded mass voltammograms different oxidation products of the model compound metoprolol were identified via their accurate masses, *e.g.* *N*-dealkylated or hydroxylated products. Especially hydroxylations are typically possible at different molecular positions. The mobility voltammograms for each *m/z* reveal isomer formation regarding hydroxylation of metoprolol. Two different hydroxylation products can be distinguished in the mobility voltammograms via differences regarding their reduced mobility referring to different collisional cross sections. Additionally MS/MS experiments after TIMS separation were carried out in order to give structural information of both compounds. Further fragmentation experiments after HPLC separation show the same results in fragmentation pattern and isomer distinction. Accordingly, by coupling EC with TIMS the analysis time can be reduced by at least the factor of five compared to conventional HPLC separations.

Conclusions

The ability to distinguish between different electrochemical metabolite isomers proves that EC/TIMS/ESI-HR-MS is a fast and powerful tool for the simulation of xenobiotic phase I metabolism including identification of the metabolites via MS/MS experiments. The time-consuming HPLC separation development process was replaced by the implementation of TIMS resulting in a faster analysis.

Novel Aspect

The new and fast separation technique TIMS has been implemented in the electrochemical simulation of phase I metabolism in order to investigate isomer formation.

ThP-149 / PRION PROTEIN CONFORMATIONAL LANDSCAPE STUDIED BY MASS SPECTROMETRY AND ION MOBILITY

Frédéric Halgand⁽¹⁾ - Guillaume Van der Rest⁽¹⁾ - Human Rezaei⁽²⁾

CNRS-Université Paris Sud, Laboratoire de Chimie Physique, Orsay, France (1) - INRA, Laboratoire VIM-MAP2, Jouy-en-Josas, France (2)

Keywords: prion, native mass spectrometry, ion mobility, conformational landscape

Introduction

MS has proved to be an efficient way to analyze “native” complexes. More recently, the commercial availability of instruments coupled to ion-mobility devices has led to renewed interest of structural biologists with access to conformational landscape of proteins. Recent studies suggest that different prion conformers might coexist in solution. We therefore decided to screen the ovine and human PrP monomers using ion mobility coupled to mass spectrometry.

Methods

Shortly, ion mobility and MS data were recorded under direct infusion of the prion protein in TEAA 20 mM pH 3.3 buffer on a Synapt G2-Si. Then collision cross section calibration was performed following the procedure of Bush et al for data comparison. For data mining we used an automated data extraction pipeline for which we developed a Python/Qt script base interface.

Results

Data showed that at least three PrP conformers are observed in the gas phase. PrP monomers are known to lead to the formation of oligomeric species in specific conditions not compatible with MS. We developed a SEC-IMS-MS setup with the aim to study the oligomers produced in these conditions. The development of this methodology allowed us to achieve CCS calibration with standards and successfully observed the small (O3) oligomeric species (O3 ~9-mer). One key observation in this process was that the relative abundance

of the gas phase monomeric conformers changed upon SEC elution. This suggests that conformer family abundances were related to protein concentration. Thus we probed the effect of protein concentration upon conformer relative abundances. Results allowed us to clearly demonstrate that monomer concentration has an effect on the ratio of conformers present in solution, which could be related to the occurrence of structural information transfer occurring through the formation of discrete polymeric objects.

Conclusions

Here we demonstrated that ion mobility coupled to mass spectrometry can be a highly efficient tool, not only to qualitatively describe the conformational landscape but also to quantitatively follow conformer dynamics. Our results opened new ways for biological interpretation of prion pathologies.

Novel Aspect

First evidence of conformers in the gas phase that keep part of solution properties and helped deciphering prion pathologies.

References

1. Van der Rest, G., F. Halgand, Size Exclusion Chromatography-Ion Mobility-Mass Spectrometry Coupling: a Step Toward Structural Biology. *J Am Soc Mass Spectrom*, 28, 2519-2522 (2017).
2. Van der Rest, G., H. Rezaei, F. Halgand, Monitoring Conformational Landscape of Ovine Prion Protein Monomer Using Ion Mobility Coupled to Mass Spectrometry. *J Am Soc Mass Spectrom*, 28, 303-314 (2017).

ThP-150 / SECONDARY STRUCTURE OF IONIZED PEPTOIDS BY ION MOBILITY MASS SPECTROMETRY

Emilie Halin - Sébastien Hoyas - Julien De Winter - Vincent Lemaire - Jérôme Cornil - Sophie Laurent - Pascal Gerbaux

UMONS, Université de Mons, Mons

Keywords: Peptoids, Secondary Structures, Ion Mobility

Introduction

Peptoids are peptidomimetic polymers of *N*-substituted glycine units that differ from amino acid building blocks by the side chain position present on the nitrogen atom [1]. Peptoids can adopt secondary structures, such as helix, which are dependent on *cis/trans* amide conformation, itself related to the composition [2]. Nowadays a great care is raised to the peptoid structure especially since the helicoidal folding possibly induces enantioselective properties [3].

Methods

(*S*)-phenylethyl peptoids are synthesized using a well-described solid-phase protocol involving a series of repeated acylation and nucleophilic displacement steps [4]. Mass spectrometry experiments have been completed using the Waters Synapt G2-Si instrument integrating the traveling wave ion mobility technology. The interpretation of experimental collisional cross sections is supported by Molecular Modeling with optimized force field.

Results

To try deciphering the gas-phase secondary structure of ionized peptoids, Ion Mobility Mass Spectrometry (IMMS) and Molecular Modeling techniques are combined. The helicoidal structure adopted by (*S*)-phenylethylamine homopeptoids (Nspe) is confirmed in solution justifying the choice of this side chain as model compound [5]. The Collisional Cross Sections (CCS) of protonated Nspe peptoids from one to fifteen

monomer units are measured. Interestingly, when the CCS values are reported in function of the polymerization degree, a linear correlation appears, potentially attributable to the helicoidal conformation. However, comparing the CCS_{exp} and the CCS_{th} and analyzing the relative energies, a charge-induced loop structure seems to be a more relevant gas phase structure. Acetylation of terminal amine function is therefore envisaged to hinder the specific charge positioning at the amine terminal function. The IMMS results obtained for AcNspe peptoids shows a totally new trend.

Conclusions

Peptoids are an emergent class of biomimetic polymers [1]. Numerous in-solution studies describe their ability to fold as a helix despite the fact that the amide conformational equilibrium drastically complicates the data interpretation [2]. Mass spectrometry represents an elegant alternative although the ionization/desolvation processes may denature the solution phase structures.

Novel Aspect

In contrast to peptides, the H-bond suppression in peptoids induced by the side chain position questions all the considerations about the structure conservation upon ESI.

References

1. Zuckermann R. N., Kerr J. M., Kent S. B. H., Moos W. H., Journal of the American Chemical Society, 114, 10646-10647 (1992).
2. Wu C. W., Kirshenbaum K., Sanborn T. J., Patch J. A., Huang K., Dill K. A., Zuckermann R. N., Barron A. E., Journal of the American Chemical Society, 125, 13525-13530 (2003).
3. Wu H., Liang T., Yin C., Jin Y., Ke Y., Liang X., The Analyst, 136, 4409–4411 (2011).
4. Tran H., Gael S. L., Connolly M. D., Zuckermann R. N., Journal of Visualized Experiments, 1-5 (2011).
5. Armand P., Kirshenbaum K., Goldsmith R. A., Farr-Jones S., Barron A. E., Truong K. T. V., Dill K. A., Mierke D. F., Cohen F. E., Zuckermann R. N., Bradley E. K., Proceedings of the National Academy of Sciences, 4309-4314 (1998).

ThP-151 / COLLISION CROSS SECTION (CCS) AS A COMPLEMENTARY PARAMETER TO CHARACTERIZE HUMAN AND VETERINARY DRUGS: CREATION OF A DATABASE

Carmen Tejada-Casado⁽¹⁾ - *Maykel Hernández-Mesa*⁽²⁾ - *Fabrice Monteau*⁽²⁾ - *Francisco J. Lara*⁽¹⁾ - *Monsalud del Olmo-Iruela*⁽¹⁾ - *Ana M. García-Campaña*⁽¹⁾ - *Bruno Le Bizec*⁽²⁾ - *Gaud Dervilly-Pinel*⁽²⁾

University of Granada, Department of Analytical Chemistry - Research group FQM302, Granada (1) - École nationale vétérinaire, agroalimentaire et de l'alimentation de Nantes-Atlantique (Oniris), Laboratoire d'Etude des Résidus et Contaminants dans les Aliments (LABERCA), Nantes (2)

Keywords: Ion mobility mass spectrometry; human and veterinary drugs; fragmentation patterns, urine samples

Introduction

Ion mobility spectrometry (IMS) in combination with mass spectrometry (MS) provides complementary information to m/z , the so-called collision-cross-section (CCS, Ω (Å²)), which may be useful for the identification of human and veterinary drugs [1]. So far, only a reduced number of these drugs has been characterized in terms of CCS, limiting its application as identification parameter [2,3]. Moreover, in the field of metabolomics where a wide number of compounds from different chemical families are detected, the creation of CCS databases is crucial for molecular identification.

Methods

IMS-MS analysis were carried out on a hybrid quadrupole/traveling wave ion mobility/orthogonal

acceleration time-of-flight (ToF) geometry instrument (Synapt G2-S HDMS). CCS characterization of human and veterinary drugs was accomplished underflow injection analysis (FIA) and positive electrospray ionization (ESI) conditions. Nitrogen was used as drift gas in the ion mobility cell. In this work, 7 families of drugs (*i.e.* 18 benzimidazoles, 11 5-nitroimidazoles, 11 aminoglycosides, 19 quinolones, 18 β -lactams, 10 sulfonamides and 5 tetracyclines) have been identified and characterized in terms of m/z and CCS.

Results

CCS is a molecular characteristic closely related to m/z , so linear correlation between both parameters is expected for compounds belonging to the same chemical family and with similar structural composition. The developed data set includes the CCS of 173 ions considering both $[M+H]^+$ and $[M+Na]^+$ species. High correlation between m/z and CCS is observed for $[M+H]^+$ ($r=0.9756$, $n=91$) and $[M+Na]^+$ ($r=0.9558$, $n=82$) ions. As expected, CCS values for sodium adducts are greater than for protonated molecules because they present higher molecular weight. However, sodium adducts of aminoglycosides, β -lactams, and of several quinolones and benzimidazoles, present lower CCS values than expected, that can be interpreted as these ions are more compact than their related protonated molecules. The main fragment ions for most compounds have also been characterized in terms of CCS and included in the database as additional interesting piece of information. Finally, 11 veterinary drugs in bovine urine samples were analyzed showing no influence of matrix on the CCS values.

Conclusions

This work proposes a CCS database for human and veterinary drugs. It includes 92 compounds of which 37 of them have been characterized in terms of CCS for the first time. In traditional LC-MS workflows, this database can be used for identification proposes together with other well established identification criteria such as precursor ion, accurate mass, fragment ion, isotopic pattern and retention time.

Novel Aspect

Additional information for the identification of a wide range of contaminants and other residues is provided as required by new trends in global and /or untargeted analysis.

References

1. Laphorn C., Pullen F., Chowdhry B. Z., Mass Spectrom. Rev. 32, 43-71 (2013).
2. Hernández-Mesa M., Escourrou A., Monteau F., Le Bizec B., Dervilly-Pinel G., TrAC Trends Anal. Chem. 94, 39-53 (2017).
3. Hernández-Mesa M., Le Bizec B., Monteau F., García-Campaña A. M., Dervilly-Pinel G., Anal. Chem. 90, 4616-4625 (2018).

ThP-152 / COLLISION CROSS SECTION AS A NOVEL PARAMETER FOR THE CHARACTERIZATION OF STEROIDS

Maykel Hernández-Mesa⁽¹⁾ - ***Bruno Le Bizec***⁽¹⁾ - ***Fabrice Monteau***⁽¹⁾ - ***Ana M. García-Campaña***⁽²⁾ - ***Gaud Dervilly-Pinel***⁽¹⁾

École nationale vétérinaire, agroalimentaire et de l'alimentation de Nantes-Atlantique (Oniris), Laboratoire d'Etudes des Résidus et Contaminants dans les Aliments (LABERCA), INRA UMR 1329, Nantes (1) - University of Granada, Department of Analytical Chemistry, Research Group FQM302, Granada (2)

Keywords: collision cross section, ion mobility spectrometry, mass spectrometry, steroids, biological samples

Introduction

Steroids play important roles in biological processes, thus, steroid profiling enhances the knowledge about metabolic mechanisms and the effects of chemical exposures on human health [1]. The integration of ion

mobility spectrometry (IM) in LC-MS workflows provides additional information to steroid profiling by adding an extra separation dimension and measuring the molecular collision cross section (CCS) [2].

Methods

A total of 300 compounds belonging to different steroid families (*i.e.* androgens, estrogens, progestogens and corticosteroids) and involving 1080 ions (*i.e.* protonated and deprotonated molecules, fragment ions, sodium adducts, dimers) have been identified and characterized in terms of m/z and CCS. A hybrid quadrupole/traveling wave ion mobility/time-of-flight mass spectrometry instrument (Synapt G2-S HDMS) was employed for CCS characterization. Urine samples containing steroid glucuronides and sulfates have been analyzed by LC-ESI-IM-MS.

Results

High correlation has been observed between m/z and CCS for steroids. However, several steroids with large moieties such as steroid glucuronides and steroid esters resulted in more compact or elongated molecules than predicted from the large set of steroids analyzed. In addition, steroids present different structural properties depending on the number of carbons of their steroid skeleton with consequences on their spatial size as evidenced by the CCSs observed. Structural differences of $22.8 \pm 5.9 \text{ \AA}^2$ have been observed between the protonated molecule and the sodium adduct of steroids. Certain protonated molecules and sodium adducts present a CCS deviated from the expected value which results in valuable information for steroid identification in untargeted analyses. Moreover, several isomeric steroid pairs (*e.g.* 5β -androstane-3,17-dione and 5α -androstane-3,17-dione) have been separated based on their CCS differences.

Conclusions

CCS provides additional information to retention time and mass spectra for supporting steroid profiling and gives more confidence to the results obtained in the routine analysis of steroids. The CCS of steroids in urine samples can be measured with high confidence considering a variation $< 2\%$ in relation to the values obtained in standards.

Novel Aspect

This work proposes the first CCS database for 300 steroids (*i.e.* endogenous, including phase I and phase II metabolites, and exogenous synthetic compounds).

References

1. Jeanneret F., Tonoli D., Rossier M. F., Saugy M., Boccard J., Rudaz S., J. Chromatogr. A 1430, 97-112 (2016).
2. Ewing M. A., Glover M. S., Clemmer D. E., J. Chromatogr. A 1439, 3-25 (2016).

ThP-153 / APPLICATION OF CO₂ DRIFT GAS IN TRAVELLING WAVE ION MOBILITY MASS SPECTROMETRY FOR THE STRUCTURAL ELUCIDATION OF ENVIRONMENTALLY RELEVANT ISOBARIC ISOMER COMPOUNDS

Yoshiteru linuma⁽¹⁾ - Eisuke Hayakawa⁽¹⁾ - Satoshi Inomata⁽²⁾ - Kei Sato⁽²⁾ - Hong Ren⁽³⁾ - Siyao Yue⁽³⁾ - Pingqing Fu⁽³⁾

Okinawa Institute of Science and Technology, Okinawa Institute of Science and Technology, Onna (1) - National Institute for Environmental Studies, National Institute for Environmental Studies, Tsukuba (2) - Tianjin University, Tianjin University, Tianjin (3)

Keywords: travelling wave ion mobility, isobaric isomers, polarizable drift gases, environmental samples

Introduction

The oxidation of volatile organic compounds in the atmosphere leads to secondary organic aerosol

compounds (SOA) that partition into ambient particulate matters. Owing to lack of authentic standards, their structural analysis mostly relies on MS_n experiments though the information tends to be limited to the presence of functional groups. Here we show the application of travelling wave mass spectrometry with polarizable drift gas for their structural analysis.

Methods

The extracts of PM samples collected in Chengdu, China were analysed using a Waters Acquity UPLC coupled to SYNAPT G2 HDMS (UPLC/(-)ESI-IMS-QTOFMS). An Acquity UPLC HSS T3 column (100 x 2.1 mm, 1.8 micrometer) was used for the separation of the sample. The ability to separate small isobaric isomer compounds was compared for N₂ drift gas and CO₂ drift gas.

Results

Among SOA compounds detected in the chromatograms, we focus on a series of isobaric isomer compounds detected at m/z 294 (m/z 294.0647, C₁₀H₁₆NO₇⁻). Based on the number of carbon and hydrogen atoms, they are thought to originate from monoterpenes. Common m/z values in MS₂ spectra were m/z 247 (C₁₀H₁₅O₅S⁻), and m/z 95 (SO₄⁻). Two of the isomers showed m/z 231 (C₁₀H₁₅O₄S⁻), whereas the one peak showed m/z 233 (C₈H₁₁NO₅⁻) instead. The presence of m/z 95 (SO₄⁻) in MS₂ suggests that precursor monoterpenes contain a tertiary carbon where a sulfate group can be attached. In the drift time of CO₂-IMS-MS, the m/z 247 fragment in the peak at 9 min showed a larger drift time (6.8 msec) than the peak eluting at 8.5 min (6.5 msec). This difference was not visible for the N₂ drift gas. This indicates that the CCS of the precursor monoterpene is larger for this peak than the other peak. While this information alone is not sufficient to determine their structures or precursors, it provides an important piece of information to narrow down their possible structures.

Conclusions

In this work, we showed that IMS-MS separation with CO₂ drift gas can be useful obtaining additional structural information to MS₂ data. The technique is potentially useful in obtaining structural information for other environmental compounds that have no authentic standard compounds.

Novel Aspect

Application of highly polarizable drift gas in IMS-MS to elucidate structures of unknown small molecules.

ThP-154 / PROFILES OF ARRIVAL TIME DISTRIBUTIONS IN DIFFERENT DRIFT GASES AS A TOOL FOR IDENTIFICATION OF ISOMERS

Lucie Borovcová⁽¹⁾ - *Jana Ondruchová*⁽¹⁾ - *Vladimír Havlíček*⁽²⁾ - *Karel Lemr*⁽¹⁾

Palacký University, Regional Centre of Advanced Technologies and Materials, Department of Analytical Chemistry, Faculty of Science, 17. listopadu 12, Olomouc (1) - Institute of Microbiology ASCR, Laboratory of Molecular Structure Characterization, Vídeňská 1083, Prague (2)

Keywords: ion mobility, mass spectrometry, isomers, new psychoactive substances, saccharides

Introduction

Ion mobility coupled to mass spectrometry improves selectivity[1] but available instruments provide relatively low resolving power. Change of drift gas can enhance separation[2] or peak fitting can reveal isomers in their mixture [3]. Although strongly overlapped peaks were distinguished by the fitting very similar peak shapes can limit its efficiency. Three drift gases were used to distinguish isomers comparing not only their drift times but also peak profiles.

Methods

Experiments were carried out using a Q-TOF instrument equipped with ion mobility (Synapt G2-S, Waters,

Manchester, UK). A custom-made gas controller (Chromservis, Prague, Czech Republic) allowed switching between different drift gases (N₂, Ar or CO₂). Sample solutions were directly infused into an electrospray ion source. Isobaric or isomeric new psychoactive substances and isomeric hyaluronan-derived oligosaccharides were used as testing compounds.

Results

Using three drift gases (N₂, Ar, CO₂), ion mobilograms of standards were acquired to get the reference drift time values and arrival time distribution (ATD) profiles. In some cases the change of drift gas allowed separation of isomeric or isobaric analytes and individual peaks were observed in mobilograms. Some critical pairs were not separated sufficiently in any tested drift gas. Nevertheless, differences in ATD (peak) profiles of some unseparated analytes were evident which allowed to distinguish them. For example for drugs of abuse, two isomers (3-FMC, 4-FMC) were well separated from isobaric 2C-H in N₂ but both isomers showed practically identical drift time values as well as ATD profiles and could not be identified. Applying CO₂, the isomers were not separated sufficiently but their ATD profiles were different and helpful for the identification. ATD profiles of some individual compounds could also be extracted from a mobilogram of their mixture and peak areas corresponding to individual compounds were determined.

Conclusions

In ion mobility, individual analytes are identified using their drift times or collision cross sections. Selectivity can be changed by applying different drift gases but drift times of compounds can still remain too close. Even in such cases, ATD profiles of individual compounds can be characteristic and useful for the identification. Since ATD profiles of analytes can be similar in one drift gas application of different gases can improve the identification ability.

Novel Aspect

In addition to drift time values, ATD profiles acquired in different drift gases were successfully used to identify isomers.

References

1. Laphorn C., Pullen F., Chowdhry B. Z., Mass Spectrom. Rev., 32, 43 (2013).
2. Howdle M. D., Eckers C., Laures A. M. -F., Creaser C. S., Int. J. Mass Spectrom., 298, 72 (2010).
3. Borovcová L., Hermannová M., Pauk V., Šimek M., Havlíček V., Lemr K., Anal. Chim. Acta, 981, 71 (2017).

Acknowledgement:

Financial support from the Ministry of Education, Youth and Sports of the Czech Republic (LO1305) and Palacký University in Olomouc (IGA_PrF_2017_020, IGA_PrF_2018_027).

ThP-155 / ANALYSING DISEASE TRIGGERING AMYLOIDS: COMBINING ESI-IMS AND LILBID-MS

Tobias Lieblein - Rene Zangl - Janosch Martin - Nina Morgner

Institute for Physical and Theoretical Chemistry, Goethe-University, Frankfurt/Main

Keywords: Alzheimer's disease, beta-amyloid, structural changes by ligands, ion-mobility, LILBID-MS

Introduction

Spontaneous aggregation of proteins can be the reason for a diversity of diseases *e.g.* Alzheimer's Disease (AD) and Parkinson. The proteins involved in these diseases (beta-amyloid (A β 42) and α -Synuclein, respectively) show the ability of spontaneous aggregation, leading to fibrils formation [1,2]. It is now generally accepted that the small oligomers rather than fibrils are the toxic species in the emergence of those diseases [3].

Methods

Laser-induced-liquid-bead-ion-desorption (LILBID) MS is an ionization technique which has the ability to conserve non-covalent interactions. Thereby, tiny droplets are transferred into vacuum and irradiated by IR laser pulses. This leads to release of the analyte-ions which are then mass-analysed by time-of-flight. Electrospray-ionization (ESI) combined with IMS allows to correlate the oligomerization of A β 42 to structural changes of the amyloid peptide.

Results

We present the differences, advantages and disadvantages of both methods for studying amyloids. By combining the two complementary techniques it is possible to analyze oligomerization, ligand binding and structural changes of the aggregating proteins in dependence of different ligands for a better understanding of the amyloid's behavior. We use LILBID to study the kinetics of the oligomerization of A β 42 and α -Synuclein by time-resolved measurements and the effect of molecules which influence the oligomerization. We show examples of molecules which inhibit or even reverse the oligomerization. Furthermore, we use ESI MS and IMS to study the peptide structure. We are able to correlate experimentally determined CCS values to the structure of aggregating peptides to detect changes of the structure induced by ligand molecules. Thereby we can gain insights into amyloid formation as well as the effect of different kinds of ligands.

Conclusions

LILBID-MS is a well-suited tool to study the kinetics of amyloid behavior allowing to follow the oligomerization of amyloids in terms of disease-relevant oligomers in dependence of the presence of inhibitors. ESI-IMS can give complementary information on how different binding partners affect the conformation of A β 42. By correlating both methods we can gain an insight on the mechanism of oligomerization.

Novel Aspect

Combining two mass spectrometry techniques to gain orthogonal information (structure and oligomeric state) of amyloidogenic peptides

References

1. Pike, C. J. ; Burdick, D. ; Walencewicz, a J. ; Glabe, C. G. ; Cotman, C. W. J. Neurosci. 13, 1676 (1993).
2. Lee, V. M. ; Trojanowski, J. Q. Neuron. 52, 33 (2006).
3. Hamley, I. W. Chem. Rev. 112, 5147 (2012).

ThP-156 / DIFFERENTIAL MOBILITY SEPARATION (DMS)-BASED SEPARATION OF BILE ACID ISOMERS

Dietrich Merkel⁽¹⁾ - Cyrus Papan⁽¹⁾ - Volkhard Kaefer⁽²⁾ - Heike Bähre⁽²⁾ - Jörg Dojahn⁽¹⁾

SCIEX Germany, Application Support, Darmstadt (1) - Medizinische Hochschule Hannover, Zentrale Forschungseinrichtung Metabolomics, Hannover (2)

Keywords: Differential Mobility Separation Bile Acids

Introduction

Bile acids are involved in a wide range of biological functions including lipid resorption, immunological functions and metabolic regulation. Through metabolic transformations, isomeric and isobaric variants are generated, which makes the unequivocal identification and quantification of individual chemical species difficult.

Here we introduce Differential Mobility Separation (DMS) as a methodology for the separation of bile-acid isomers. DMS is an ion mobility technology which separates molecules based on their dipolar moment. We show DMS separation in conjunction with chromatographic separation (LC-DMS-MS) as well as with direct infusion (DMS-MS). While the combination with chromatographic separation may improve selectivity, the

separation power of SelexION is sufficient for a clear separation of isomers, allowing for infusion-based fast quantification without the need for LC development.

Methods

Bile acid standards were prepared in 10% DMSO at a stock solution of 1 mM. Two groups were used: taurodeoxycholic acid (TDCA), taurochenodeoxycholic acid (TCDCa), tauroursodeoxycholic acid (TUDCA) with monoisotopic mass of (499.2967) and formula ($C_{26}H_{45}NO_6S$) and glycodeoxycholic acid (GDCA), glycochenodeoxycholic acid (GCDCA) and glyoursodeoxycholic acid (GUDCA) with monoisotopic mass of (449.3141) and formula ($C_{26}H_{43}NO_5$). For direct infusion, compounds were prepared in water/methanol 50/50, 0.1% FA at a concentration of 1 μ M. For LC-separation preparation was done in the same solvent at a concentration of 0.1 nmol/ml with 1 μ l injection volume. The samples were measured on a QTRAP[®] 6500 LC-MSMS system equipped with a SelexION[®] (DMS) device and coupled to an ExionLC[™] system for the LC experiments.

Preliminary Results

Using the DMS technology, we directly infused the bile acid standards individually and as mixtures to determine the CoV values of the different isomers. To this end, the compensation voltage was ramped over a range between -30V to 0V. The results showed almost baseline separation of the different isomers.

The next aim was to combine DMS with chromatographic separation. The combination of DMS with chromatography has several advantages: Firstly, DMS showed a marked reduction of chemical background for better quantification, resulting in improvement of the signal to noise ratio (give a value NNN). A further advantage is less requirement for development of chromatographic separation, which was shown by the elimination of isomers from extracted ion chromatograms. This also eliminated the requirement for retention time assignment and less reliance on compound specific fragments which may be low abundant. For example, DMS allows to use a highly abundant fragment without the need for compound specificity.

Novel Aspect

Differential mobility separation, increased selectivity of bile acid separation by LC-DMS-MS

ThP-157 / ADDITION OF COLLISION CROSS SECTIONS TO SEARCHABLE LIBRARIES: A ROBUST APPROACH TO ENHANCE CUMULATIVE SPECIFICITY IN SCREENING ASSAYS

Russell Mortishire-Smith⁽¹⁾ - ***Mike McCullagh***⁽¹⁾ - ***Mélanie Herman***⁽²⁾ - ***Els Van Hoeck***⁽²⁾ - ***Séverine Gosciny***⁽²⁾

Waters, Scientific Operations, Wilmslow (1) - Sciensano, Department of Food, Medicines and Consumer Safety, Brussels (2)

Keywords: Collision Cross Section, Ion Mobility, Steviol Glycosides, Screening Libraries, CCS Fingerprints

Introduction

CCS data bases are being developed across multiple application areas (1-6). A method using Ion Mobility Mass Spectrometry has been developed to screen/distinguish steviol glycoside isomers. A quantitative assessment of long-term CCS measurement stability and an inter site comparison is presented. In addition to illustrating how CCS fingerprints can be produced to enhance cumulative identification specificity within screening libraries.

Methods

Reversed phase separation was performed (2 μ l/min), using a microfluidic PCA (post column addition) separation device (40°C), incorporating a 1.7 μ m, ACQUITY UPLC BEH C18 stationary phase in a 150 μ m diameter separation channel, with integrated ESI emitter (-2.6 kV). The source enclosure was interfaced to a SYNAPT[®] G2-Si mass spectrometer and generic ion mobility separation performed in negative ion HDMSE mode, at 20000 FWHM (m/z 50-2000).

Results

55 food commodity samples labeled as containing stevioside sweeteners were screened using a steviol glycoside CCS library. Analyses produced CCS values, enabling the unequivocal identification of steviol glycosides and isomeric pairs. Retrospective data analysis revealed coeluting isomeric species comprising, (steviol glycosides, formate adducts, labile fragment ions, doubly de-protonated dimers, multiply charged species) and have been differentiated using CCS, generating specific isomeric analyte fingerprints.

Short term precision/repeatability CCS measurement RSD's between 0.11 and 0.21% were obtained. For an inter-site comparison (UK/Belgium), the CCS difference was <0.7%. For long term performance, over a 45 days (24/7, one calibration), samples comprised of solvent standards (1-1000 pg/ μ L) and a commodity subset. 18,974 steviol glycoside detections were measured. The determined CCS error was <1%, with RMS error of 0.26%.

Conclusions

Inter-site and short/long term precision of CCS measurements indicate that CCS is a reliable metric that can be added to searchable screening libraries, enabling an additional complimentary metric to be used in conjunction with retention time and accurate mass measurement. CCS fingerprints can be used as an additional metric to enhance identification specificity, even where no product ion information is available and retention time shifts occur.

Novel Aspect

Micro UHPLC MS sensitivity and CCS retrospective data review, enabled multiple analyte specific isomeric species (CCS fingerprints) to enhance searchable databases.

References

1. Pulkrabová J., Tomaniová M., (6th International Symposium on Recent Advances in Food Analysis (RAFA), Eds. (2013). November 5–8, Prague, Czech Republic.
2. Gallagher R., Pattison C., Pickup K., Chipperfield J., McCullagh M., Douce D., Proc. 62st ASMS Conf. Mass Spectrometry and Allied Topics, Baltimore, USA, June (2014).
3. Righetti L., Bergmann A., Galaverna G., Rolfsson O., Paglia G., Dall'Asta C., Analytica Chimica Acta (2018) (article in press).
4. Lipok C., Jörg Hippler L., Schmitz O. J., Journal of Chromatography A, 1536, 50–57 (2018).
5. Bauer A., Kuballa J., Rohn S., Jantzen E., Luetjohann J., Journal of Separation Science doi:10.1002/jssc.201701059. (article in press).
6. Hernández-Mesa M., Le Bizec B., Monteau F., García-Campaña A. M., Dervilly-Pinel G., Anal. Chem. 90, 4616-4625 (2018).

ThP-158 / THE USE OF CCS-ENABLED LIBRARIES FOR DRUG METABOLISM APPLICATIONS

Russell Mortishire-Smith⁽¹⁾ - ***Daniel Veyel***⁽²⁾ - ***Bartłomiej Krawczyk***⁽²⁾ - ***Hans Vissers***⁽¹⁾

Waters, Applied Informatics, Wilmslow, United Kingdom (1) - Boehringer Ingelheim, Drug Discovery Sciences, Ingelheim am Rhein, Germany (2)

Keywords: CCS, Libraries, DMPK, Unknowns, Matrices

Introduction

A primary challenge in characterizing metabolic fate is distinguishing drug-related material in the presence of the complex matrix. A commonly used approach is to filter components based on parent drug properties. However, identification of drug-related, unexpected metabolites remains a limiting issue. By characterizing the native complement of biological matrices in a library based on m/z , CCS and product ions, we can flag components in samples that are unlikely to be drug metabolites.

Methods

Hepatocyte incubations were analysed using ion mobility enabled data independent acquisition (HDMSE) LC/MS data were acquired with two different reversed phase gradient methods. The data were co-detected across all samples in the m/z , RT, and CCS dimensions. An aggregate peak list file was developed from which an initial query library was created. Through an iterative search and (Kendrick) filtering process, the content of the known-unknown library was curated and appropriate query parameters derived.

Results

Five technical replicate LC-MS experiments of each incubation were conducted per gradient method. Four of the five replicates were used for library creation and the remaining replicate was used for testing and evaluation. The initial library comprised components that were consistently detected across all runs with regard to m/z , retention time and drift time reproducibility. No restrictions were applied at this stage with respect to intensity and/or number of fragment ions per product ion spectrum. Appropriate search tolerances were readily obtained by considering every possible match parameters, resulting in a set of parameters that are equivalent to those typically applied in studies where the structure and/or elemental composition of the target molecules are known. The curation process also allowed for the estimation of the contribution of each analytical parameter, either individually or combined, to the overall specificity of a compound search. In addition, two methods were considered and evaluated that express normalized specificity metrics either as a function of the number of detections or curated library entries, respectively.

Conclusions

The application of the known-unknown library and the ability to estimate FDR will be demonstrated by hepatocytes incubation matrices spiked with known compounds and metabolites. The development and evaluation additional known-unknown libraries of other biological matrices is currently underway as well as the identification of commonly detected unknown compounds in different biological drug incubation matrices.

Novel Aspect

Our approach shows the benefit of integrating CCS values even for unknown compounds in known-unknown libraries which has the potential to accelerate metabolite identification in DMPK workflows.

ThP-159 / SEPARATION OF INTACT PROTEIN ISOFORMS WITH DIFFERENTIAL MOBILITY

Kerstin Pohl⁽¹⁾ - *Chris Nortcliffe*⁽²⁾ - *Ferran Sanchez*⁽³⁾ - *Jörg Dojahn*⁽¹⁾ - *Annu Uppal*⁽⁴⁾ - *Sibylle Heidelberger*⁽²⁾ - *Yves J. C. Le Blanc*⁽⁵⁾

SCIEX, Biologics, Darmstadt (1) - SCIEX, Biologics, Warrington (2) - SCIEX, Biologics, Barcelona (3) - SCIEX, Biologics, Gurugram (4) - SCIEX, R&D, Concord (5)

Keywords: Protein Structure, Ion Mobility, HDX, DMS, Gas Phase Structure

Introduction

Native electrospray offers the ability to maintain structural characteristic of the proteins prior to MS analysis. In recent years, the coupling of conventional mobility has also proven to be useful at interrogating the various conformation that may exist in the gas phase associated with intact proteins. More recently, differential mobility separation (DMS) has also been used to effectively separate isoforms and isomers prior to analysis. Here we report of the use of DMS to augment the separation of intact protein under native- LC conditions.

Methods

All proteins and peptide standards were acquired from SIGMA and used without any further purification. LC was performed using a Shimadzu Nexera UPLC system operated at a flow rate of 400 μ L/min with a

2x150mm Yara 3u SEC-3000 (Phenomenex) over 15 minutes. All samples were analysed using a TripleTOF® 5600+ system equipped with a SelexION® (DMS system) using a Turbo V source.

Results

Cytochrome-c was used as model protein to assess the DMS and source condition required for successful detection of various conformations. Under 'native like' infusion conditions, it was found that the DMS can be operated at medium temperature (225°C) and varying SV values to reveal multiple transmission optima (n=3) for the +8 charge state. Prior publications have observed and attributed multiple conformers to this particular charge state. These conformers were also observed under LC conditions as long as the source temperature was kept below 300°C. Though CCS measurement may not be directly obtained from the CoV values, HDX in the gas phase after the DMS separation confirmed that the conformations had differences in exchange rate that are consistent with compact and extended structure. When the spray and solution conditions were modified to induce conformational changes in the solution phase, the CoV optima for that charge state changed. These results are consistent with conformational changes that can be detected by DMS separation.

Conclusions

Differential mobility provides an orthogonal separation to mass spectrometry and LC technology. The work presented here opens up an avenue of structural biology, with the interrogation of compact and extended conformations. The addition of gas phase HDX provides further evidence that the ions are entering the gas phase with a level of solution phase structure. Importantly the observed species are supported through previous publications.

Novel Aspect

'Native like' structure analysis with differential mobility and LC separation

ThP-160 / ION GUIDE WITH VERY STRONG AXIAL TRANSPORTING FIELDS

Sergey Poteshin - Anna Burykina

NRNU MEPHI, Molecular physics, Moscow

Keywords: Ion transport, ion transport, Collision Cell with axial field

Introduction

The main task of this investigation is the searching of geometrical and RF power supplying modes for RF-only stacked ring ion guide with traveling waves for providing very high axial transporting fields. It can be attractive for 1) high-speed, high-pressure collisional cell, 2) the lifting of ions from low to high pressure area for Quadrupole/IMS applications.

Methods

For a theoretical calculation of effective potentials for the stack of ring electrodes with RF supplying we can use comprehensive work of Berdnikov A. S. [1]. Also we calculate ions trajectories in gas with using SIMION 8.1 with own modified HS1 code. The modification allows to significantly accelerate the calculation of trajectories at the high pressure without loss of accuracy.

Results

As a basic version of stacked ring ion guide we chose a 4 times the electrode pitch which repeated along axis. Thus we needed only 4 RF power signal for applying to N=1. . 4 electrodes. RF signal is a combination high frequency V_{rf+} and low frequency (LF) signal providing traveling waves field. The LF voltages have an amplitude V_{tw} and are phase shifted by $\pi/2N$.

The inner diameter of the electrodes was chosen quite small $R_0 \sim 1$ mm, to produce straight axis field caused by TW potentials. The thickness of electrodes was 0.2 mm and the gap between of them is 0.2 mm. Numerical simulation shows the good ion transport efficiency at wide pressure range 0-100 Torr with different RF power supplying regime. At low pressure we have to apply higher amplitude with higher frequency V_{rf+} and for high pressure we can leave only LF supplying signals for the traveling waves. Also simulation shows excellent ability for a high speed axial ion transporting.

Conclusions

The ability of ion guide based on stacked ring with traveling waves to operate at wide pressure range 0-100 Torr with high axial transporting speed of 100 m/s (at 100 Torr) is shown. Thus our ion guide can operate as ions passing system from low to high pressure ~ 10 -100 Torr for realizing Quadrupole/IMS devices. Also it will be shown that the presented ion guide is effective as a very high speed at operating rate collision cell.

Novel Aspect

Represented ion guide can be used for 1) ions passing system from low to high pressure for Q/IMS devices or high pressure collision cell 2) very high speed operating rate collision cell.

References

1. Berdnikov A. S., Nauchnoe Priborostroenie (Scientific Instrumentation), 2011, 21, № 4

ThP-161 / THERMAL DENATURATION OF PROTEINS BY VARIABLE-TEMPERATURE ELECTROSPRAY IONIZATION PAIRED WITH ION MOBILITY SPECTROMETRY-MASS SPECTROMETRY

Shannon Raab - Tarick El-Baba - Daniel Woodall - David Clemmer

Indiana University, Department of Chemistry, Bloomington

Keywords: IMS-MS, vt-ESI, protein folding, structure

Introduction

Protein folding is difficult to study due to the lack of sensitive techniques that can detect all populated structures on a folding pathway. Ion mobility spectrometry-mass spectrometry (IMS-MS) paired with a variable-temperature electrospray ionization (vt-ESI) source has elucidated multiple folding structures based on changes in collision cross section (CCS) upon protein denaturation [1]. Here, this method is used to study protein folding and assembly.

Methods

Experiments were carried out using nested IMS-TOF instruments fitted with a custom vt-ESI source that maintains a solution temperature up to ionization and analysis. Wild-type and point-mutated chymotrypsin inhibitor 2 (CI2) were expressed in *E. coli* and purified by existing techniques. Avidin was purchased from Lee BioSolutions and all other reagents were purchased from Sigma Aldrich. Samples were made to ~ 30 μ M with varying amounts of ammonium acetate buffer.

Results

We monitored the unfolding of CI2, a single-domain protein previously reported to fold by a two-state cooperative transition [2]. Changes in the charge state distributions suggest the melting transition follows a two-state transition; our melting temperatures (at several solution pH values) are within 1% of calorimetric data. However, the mobility distributions within each charge state reveal several structures that behave differently with increased temperature (*i.e.*, multiple folding states). Point-mutated CI2 variants were melted and compared to the wild-type protein to determine the contribution of specific residues to the overall folding behavior. Preliminary experiments studying the homotetrameric protein avidin show

interesting melting behavior. For example, the compact tetramer is favored until $\sim 82^{\circ}\text{C}$; beyond this point the tetramer appears to disassociate into folded monomers. Further increases in temperature result in unfolding of the monomers by what appears as a two-state process.

Conclusions

Despite the two-state behavior observed in the mass spectral data, CCS structures provide a clear indication that two-state behavior is a composite of many structures. CCS data show a minimum of two native states and four denatured states for CI2. The results for avidin show an example of multiple cooperative transitions hidden within a complex dissociation and unfolding transition – the tetramer first dissociated, followed by monomer unfolding.

Novel Aspect

Multiple structures are involved in a protein typically thought to unfold via a two-state transition.

References

1. El-Baba, T. J., Woodall, D. W., Raab, S. A., Fuller, D. R. Laganowsky, A., Russell, D. H., Clemmer, D. E., Journal of the American Chemical Society, 139, 6306-6309 (2017).
2. Jackson, S. E., Fersht, A. R., Biochemistry, 30, 10428-10435 (1991).

ThP-162 / SEPARATION OF ISOMERS WITH ULTRA-HIGH-RESOLUTION ION MOBILITY SPECTROMETRY

Christian-Robert Raddatz - Ansgar T. Kirk - Stefan Zimmermann

Leibniz Universität Hannover, Institute of Electrical Engineering and Measurement Technology, Sensors and Measurement Technology, Hannover

Keywords: ultra-high resolution drift tube ion mobility spectrometry, isotopes, isotopologues, collision cross section, atmospheric pressure ionization

Introduction

In drift tube ion mobility spectrometry (IMS), ions are driven through a drift tube by a homogenous electric field and collide with neutral drift gas molecules. The separation mainly depends on the ion-molecule collision cross section (CCS). Therefore, in more and more applications IMS are used for pre-separation of isomers that cannot be distinguished by their mass-to-charge ratio (m/z) in mass spectrometers (MS). One major benefit of IMS is the fast separation within a few milliseconds revealing the CCS.

Methods

The ultra-high resolution drift tube IMS used in this work has a resolving power of $R_p=250$, which is defined as the full width at half maximum (FWHM) divided by the drift time. To achieve such high resolving power in IMS a drift voltage of 25 kV is used. Nevertheless, the length of the drift tube is just 15 cm. For ionization, atmospheric pressure chemical ionization with a miniature soft X-ray source is used as well as atmospheric pressure photo ionization (APPI) with 10.6 eV.

Results

In previous work, we demonstrated the separation of several isotopologues with the instrument used. [1] Here, isotopes with small differences in CCS are separated. The percentage difference in CCS is defined as the difference of the CCS divided by the average CCS and multiplied by one hundred. Calculating the percentage difference that can be resolved with $R_p=250$ yields a minimum percentage difference in CCS of 0.8%. [2] For example, the isomers 2-picoline and 3-picoline are clearly separated using an APCI ionization source. Whereas these substances can hardly be resolved with GC-MS. [3]

Conclusions

The separation of isomeric compounds is of major interest, especially in pharmaceutical science. With ultra-high resolution drift tube IMS and IMS-MS small percentage differences in CCS can be resolved.

Novel Aspect

Separation of isotopes with differences of only 0.8 % in CCS using ultra-high-resolution DT-IMS.

References

1. A. T. Kirk, C. -R. Raddatz, S. Zimmermann, *Anal. Chem.*, 89, 1509–1515 (2017).
2. J. N. Dodds, J. C. May, J. A. McLean, *Anal. Chem.*, (2017).
3. N. P. Kulshreshtha, S. C. Moldoveanu, *J. Chrom. A*, 985, 303–312 (2003).

Acknowledgements

This work is funded by the Deutsche Forschungsgemeinschaft (DFG, German Research Foundation) – ZI 1288/4-1.

ThP-163 / DECONVOLUTION OF DIRECT INFUSION TIMS-MS ANALYSES

Peter Sander - Sven Meyer

Bruker Daltonics, R&D, Bremen

Keywords: Trapped Ion Mobility Spectrometry; deconvolution; isobaric compounds; CCS determination

Introduction

The combination of ion mobility separation with mass spectrometry offers increased capabilities for determining the structure of molecules. In these scenarios, ion mobility initially delivers a separation based on molecule sizes. Yet, often even simple molecules or molecule mixtures show ambiguous mobilities which cannot easily be resolved further, *e.g.*, by chromatography.

An algorithm has been developed to deconvolute these overlapping signals in the mobility vs. m/z domain of a direct infusion experiment. The result are compounds comprising each of a clean mass spectrum and a mobilogram. Optionally, a fragment mass spectrum is extracted.

Methods

A mixture of isobaric lipids has been analyzed using direct infusion trapped ion mobility spectrometry (TIMS) coupled with an orthogonal time-of-flight mass spectrometer.

The algorithm developed here is evaluating such analyses. It is a further development of an algorithm initially developed for deconvolution of compounds of a chromatographically separated LC/MS analysis. [1] The basic idea of this algorithm is the comparison of all extracted ion mobilograms to detect and deconvolute compounds. The key to success is the usage of fuzzy logic for the decision if any of two ions belong to the same compound or to different, partly overlapping compounds. Each compound is then associated with a clean compound mass spectrum, deconvoluted by any overlapping compounds and background subtracted. In extension to that fragment spectra at the same mobility are also associated to that compound.

Finally, the collision cross section of each compound is determined from the m/z ratio of the most intense ion and its charge state.

Results

The mobility- m/z heatmap of a measurement of a mix of lipids reveals several partly overlapping signals. The direct infusion analysis was acquired over one minute and was performed on the Bruker timsTOF Pro instrument.

Applying the algorithm delivers compounds, each of a clean mass spectrum and a mobilogram. Signals are clearly separated and contain no signals from overlapping compounds or any background signals. Since charge state and molecular weight are also determined, collision cross section (CCS) values are calculated as well and the mobilograms can be directly depicted on a CCS scale.

Summarizing the results shows that also isomeric lipids were separated, each measured with a distinct CCS value.

Conclusions

Ion mobility separation has already proved to be a valuable combination with mass spectrometry. In the task of structure elucidation however, even simple mixtures of isobaric compounds can show complex and ambiguous signals in a direct infusion experiment. This algorithm (DissectTIMS) deconvolutes these signals into compounds, each of a clean mass spectrum, a mobilogram and a CCS value.

Novel Aspect

Deconvolution of overlapping signals in a direct infusion TIMS-MS experiment delivering clean mass spectra, mobilograms and CCS values.

References

1. P. Sander, C. Baessmann, B. Schneider, G. Zurek, D. Wunderlich., Proceedings 52nd ASMS Conference (2004), Deconvolution of Complex LC/MSn Analyses

ThP-164 / ION MOBILITY MASS SPECTROMETRY FOR PROFILING AND STRUCTURAL CHARACTERIZATION OF HUMAN HIPPOCAMPUS GANGLIOSIDES

Mirela Sarbu⁽¹⁾ - *Željka Vukelić*⁽²⁾ - *David Clemmer*⁽³⁾ - *Alina Diana Zamfir*⁽¹⁾

National Institute for Research and Development in Electrochemistry and Condensed Matter, Department of Condensed Matter, Timisoara, Romania (1) - University of Zagreb Medical School, Department of Chemistry and Biochemistry, Zagreb (2) - Indiana University, Department of Chemistry, Bloomington (3)

Keywords: ion mobility separation mass spectrometry; gangliosides; human fetal hippocampus; collision-induced dissociation; cholinergic activity.

Introduction

The structural analysis of gangliosides (GGs), a particular class of glycosphingolipids highly expressed in central nervous system, experienced a new level due to the rapid advances in ion mobility separation mass spectrometry (IMS MS) [1]. In view of the direct correlation of GGs with neurogenesis and synaptogenesis, in this work we have implemented IMS MSn for the first time in the study of a highly complex native mixture originating from human fetal hippocampus [2].

Methods

The GG mixture was purified from a 17-th gestational week fetal hippocampus (FH17) using the method of Svennerholm and Fredman and modified by Vukelić et al [4]. The sample dissolved in methanol to 5 pmol/μL concentration was infused at 2μL/min flow rate and 2 kV ESI voltage into the Synapt G2s. To enhance the separation IMS wave velocity was set at 650 m/s, IMS wave height at 40 V and IMS gas flow at 90 L/min. CID was performed in the transfer cell, using energies between 10-45 eV.

Results

The ionization and detection of low abundant species was possible due to IMS that allowed the differentiation of GGs based on the carbohydrate chain length and degree of sialylation. These findings provided a more reliable insight into the expression of polysialylated structures, which might constitute

markers of brain development. By IMS MS, after six minutes of signal acquisition, no less than 134 distinct GGs were identified, a number nine times higher than the number of GGs previously reported in FH following chip nanoESI QTOF MS analysis without prior separation (only 14 structures) [3]. Moreover, IMS MS evidenced no less than 54 different species modified by either carbohydrate or non-carbohydrate type of attachments. Of these, 25 ions correspond to GG structures modified by CH₃COO⁻ attachment, novel species never reported before. CID MSⁿ experiments conducted after mobility separation provided a detailed structural characterization of CH₃COO⁻ GT1b(d18:1/18:0), a biologically relevant compound.

Conclusions

The screening of FH17 by IMS MS provided data related to the expression and role of GGs in fetal hippocampus, with a particular emphasis on their cholinergic activity by the discovery for the first time of novel species modified by CH₃COO⁻ attachment. In view of the present findings, we may conclude that IMS MS and CID MSⁿ proved their high efficiency for unequivocal detection and characterization of these low abundant species, with a significant biological role.

Novel Aspect

First implementation of IMS MS for GG determination in FH and discovery of CH₃COO⁻ modified species, correlated with the cholinergic activity of GGs in central nervous system.

References

1. Sarbu M., Robu A. C., Ghiulai R. M., Vukelić Ž., Clemmer D. E., Zamfir A. D., *Anal. Chem.*, 88, 5166-5178 (2016).
2. Sarbu M., Vukelić Ž., Clemmer D. E., Zamfir A. D., *Biochimie*, 139, 81-94 (2017).
3. Vukelić Ž., Zarei M., Peter-Katalinić J., Zamfir A. D., *J. Chromatogr. A*, 1130, 238-245 (2006).
4. Vukelić Ž., Metelmann W., Müthing J., Kos M., Peter-Katalinić J., *Biol. Chem.*, 382, 259–274 (2001).

ThP-165 / A COMPUTER CODE FOR OPTIMIZATION OF HADAMARD TRANSFORM AND FOURIER TRANSFORM ION MOBILITY SPECTROMETRY EXPERIMENTAL CONDITIONS

Anastasia Sarycheva - Alexey Adamov - Sergey Poteshin - Alexey Sysoev

National Research Nuclear University MEPhI, Molecular Physics Department, Moscow

Keywords: Hadamard transform, Fourier transform, signal-to-noise ratio, ion mobility spectrometry, signal multiplexing

Introduction

Hadamard transform was effectively applied to optical spectroscopy [1] and later was implemented in transient physical methods such as time-of-flight mass spectrometry [2] and ion mobility spectrometry [3]. Recently the method was significantly improved by pseudorandom sequence modification proposed by Belov with co-authors [4]. The modifications ensure considerable signal-to-noise ratio increase while keeping sufficiently high resolving power.

Methods

The mentioned above technique was implemented and Matlab code was developed. The code considers properties of model compounds; produces pseudorandom sequence; simulates the multiplexed spectra formation taking into account ion beam modulation, ion diffusion, space charge effects; it features ion mobility signal reconstruction. We used the code in order to consider different algorithms of ion beam modulation and ion mobility separation optimization parameters.

Results

First, we tested the code using the Fourier transform and Hadamard transform ion mobility spectrometry

data collection and treatment algorithms. The modeled multiplexed signal includes statistical variation as well as an electronic noise. The reconstructed signals for Hadamard transform demonstrate considerable rise of signal-to-noise ratio as expected. No clear advantage for signals reconstructed by Fourier transform ion mobility spectrometry data was found.

Then we used the code for optimization of pseudorandom sequencing considering signal-to-noise ratio an optimization parameter. A homologous series of tetraalkylammonium halides was chosen as the model compounds with well-defined properties in model experiments. Different multiplexing/reconstructing approaches were studied. The presentation describes main algorithms and demonstrates examples of signal-to-noise ratio improvement.

Conclusions

The proposed approach and the developed code allowed one to optimize experimental conditions in order to significantly improve signal-to-noise ratio.

Novel Aspect

Novelty of the approach is based on the algorithms allowing one to optimize multiplexing sequencing and reconstructing algorithms.

References

1. Harwit M., Sloane N. J. A., Hadamard Transform Optics, Academic Press: New York (1979).
2. Brock A., Rodrigues N., Zare R. N., Anal. Chem., 70, 3735 (1998).
3. Szumilas A. W., Ray S. J., Hieftje G. M., Anal. Chem., 78, 4474 (2006).
4. Clowers B. H., Belov M. E., Prior D. C., Danielson W. F., III, Ibrahim Y., Smith R. D., Anal. Chem., 80, 2464 (2008).

ThP-166 / FAIMS-MS-IR SPECTROSCOPY WORKFLOW: A MULTIDIMENSIONAL PLATFORM FOR THE ANALYSIS OF MOLECULAR ISOFORMS

Baptiste Schindler - Isabelle Compagnon - Franck Lepine

Institut Lumière Matière, Université de Lyon, Lyon

Keywords: FAIMS, IRMPD spectroscopy, isomers, hyphenated techniques, multidimensional analysis

Introduction

The analysis of isoforms is challenging by mass spectrometry and ion mobility spectrometry is gaining momentum for this purpose. However IMS technologies suffer from the absence of detailed structural information. On the other hand, IRMPD spectroscopy (Infrared Multiple Photon Dissociation) is able to resolve subtle isomerisms on pure sample. We propose to combine FAIMS separation, mass spectrometry and IRMPD spectroscopy for multidimensional analysis. [1]

Methods

This new instrument consists of an ultra FAIMS (Owlstone) device interfaced to a linear ion trap (LTQ XL Thermo Scientific) which was modified for IRMPD spectroscopy. The combination of FAIMS separation, mass spectrometry and infrared spectroscopy in a single instrument allows various possibilities of multidimensional analysis, especially for mixture of isomers.

Results

For a given mass-selected ion, two different modes using FAIMS separation and IRMPD spectroscopy are enabled to characterize the composition of a mixture of isomers. They are demonstrated on an isomeric mixture of paracetamol and 2-phenylglycine. In the first mode, a FAIMS (high-Field Asymmetric waveform

Ion Mobility Spectrometry) separation of the isomers is performed with a static compensation field for mass- and isomer- selective ion spectroscopy. In the second mode, the compensation field is scanned while the ions are irradiated at a fixed wavenumber allowing selective photofragmentation of a chosen species. the FAIMS-MS-IR workflow allows separation and identification of mixtures of isomers using an appropriate database of IRMPD spectra and without the use of standards for calibration purpose. Alternatively for such isomeric mixture, IRMPD spectroscopy alone yields a convoluted spectrum. By adding a FAIMS dimension, it is possible to retrieve individual spectroscopic signature for each isomer present in the mixture.

Conclusions

This configuration of IMS-MS-IR workflow constitutes an attractive alternative to more elaborate approaches for ion spectroscopy of mass-selected and isomer-selected ions. In the field of analytical chemistry, we hope that this type of multidimensional analysis will be valuable tool for broad range of biological/medical/pharmaceutical applications where the identification of isoforms is critical.

Novel Aspect

A new multidimensional analysis workflow combining FAIMS separation and IRMPD spectroscopy for analytical purpose or structural analysis.

References

1. Schindler B., Depland A. D., Renois-Predelus G., Karras G., Concina B., Celep G., Maurelli J., Loriot V., Constant E., Bredy R. et al. International Journal for Ion Mobility Spectrometry 20, 119–124 (2017).

ThP-167 / PROTOMER SEPARATION, CHARACTERIZATION AND ION-MOLECULE REACTIONS USING DIFFERENTIAL ION MOBILITY

Pascal Schorr - Dietrich A. Volmer

Humboldt University of Berlin, Analytical and environmental Chemistry, Berlin

Keywords: DMS, Protomers, Ion-molecule reactions

Introduction

Solution basicity often differs from gas phase proton affinities, generating protomers, which often differ in their CID behavior [1-3]. Information on the protomer distribution has been obtained from DMS separations [1,2]. While ion-molecule reactions have been studied in rf-only collision cells and ion traps [4], they have not been systematically investigated in DMS devices. This presentation describes controlled ion molecule reactions in a DMS cell.

Methods

Experiments were performed on an AB Sciex QTRAP 5500 system equipped with a Selexion differential ion mobility cell. The mass spectrometer was operated in positive electrospray ionization mode and samples of 4-quinolone antibiotics in methanol/water +0.1% formic acid were injected continuously at 7 $\mu\text{L}/\text{min}$ using a syringe pump.

Results

DFT calculations predict two protomers of ciprofloxacin during the ESI process: protonation on the piperazinyl moiety, specifically at N-22, and protonation at the O-atom of the C-4 keto function. These theoretically calculations were supported by DMS-MS/MS experiments. Separations of 4-quinolone protomers were achieved using DMS with a variety of chemical modifiers and were confirmed by different decay mechanisms after CID. [1]A class-selective reaction for 4-quinolone antibiotics was utilized to demonstrate the ability of performing chemical reactions in a differential ion mobility cell. Seven auxiliary

compounds as chemical reagents in the DMS cell and four 4-quinolone antibiotics were used in this study. Dehydrated reactive precursor ions were formed and annulation reactions with the reagents occurred during the residence time within the DMS. Collision energy differences of adducts and newly formed compounds clearly demonstrated formation of covalent bonds.

Conclusions

In this study we combined DFT calculations, CID mass spectra and DMS separations to describe controlled chemical reactions in a DMS device.

Novel Aspect

We report for the first time deliberate ion-molecule reactions in a DMS cell.

References

1. Kovačević B., Schorr P., Qi Y., Volmer D. A., J. Am. Soc. Mass, 25, 1974-1986 (2014).
2. Lalli P. M., Iglesias B. A., Toma H. E., de Sa G. F., Daroda R. J., Silva Filho J. C., Szulejko J. E., Araki K., N. Eberlin M. N., J. Mass. Spectrom., 47, 712 – 719 (2012).
3. Kass S. R., Tian Z., Angew. Chem. Int. Ed., 48, 1321 – 1323 (2009).
4. Colorado A., Brodbelt J., Anal. Chem., 66, 2330-2335 (1994).

ThP-168 / PROTEIN ANALYSIS USING ION MOBILITY MS AND IN-SOURCE ION ACTIVATION COMBINED WITH ECD ION FRAGMENTATION

Ruwan Kurulugama⁽¹⁾ - *Valery Voinov*⁽²⁾ - *Kenneth Newton*⁽¹⁾ - *Christian Klein*⁽¹⁾ - *Nathan Lopez*⁽²⁾ - *George Stafford*⁽¹⁾ - *Joseph Beckman*⁽²⁾ - *John Fjeldsted*⁽¹⁾

Agilent Technologies, R&D, Santa Clara, CA, United States (1) - e-MSion Inc., Corvallis, OR, United States (2)

Keywords: In-source, Protein, Activation, Mobility, ECD

Introduction

Tandem mass spectrometry is accepted as a powerful analytical technique and when used with electron capture dissociation (ECD) makes it well suited for protein identification. In-source ion activation can be used to facilitate collision induced unfolding (CIU) of proteins and when coupled to IM-MS offers information about ion energetics and structures [1]. Combining the above technologies may provide top-down protein analysis with improved sequence coverage.

Methods

An Agilent 6560 IMQTOF instrument was modified to include ion optic elements at the exit of the ion transfer capillary to improve in-source ion activation prior to ion mobility analysis. The instrument was also equipped with a prototype e-MSion ECD cell mounted at the exit of the quadrupole mass filter. Proteins were prepared under native and denatured conditions and data were analyzed using IM-MS browser, CIU suite programs, ProSight, and LcMsSpectator.

Results

Native ESI conditions for carbonic anhydrase and BSA show low charge state precursor ions ranging from +8 to +11 and +15 to +19 respectively. These native low charge state ions have folded compact structures shown by their measured small CCS values, but they can be unfolded by using in-source ion activation. These proteins under denatured ESI conditions provide higher precursor ion charge states, +12 to +38 and +22 to +58 respectively, with correspondingly elongated ion structures and larger measured CCS values. The ECD cell provided nearly complete sequence coverage of denatured proteins as large as carbonic anhydrase with preservation of post-translational modifications. Sequence coverage results will be presented for the above proteins under native and denatured ionization conditions, coupled with CIU to

facilitate more extensive fragmentation, and combined with ECD technology. Experimental data from larger proteins like monoclonal antibodies will also be analyzed and presented.

Conclusions

Electron based ion fragmentation methods, such as ECD, produce extensive ion fragmentation for readily exposed regions of proteins. The fragmentation efficiency for a protein can be enhanced by either using denatured proteins or by unfolding native proteins in gas phase prior to ECD fragmentation. This study focuses on the use of CIU to improve ECD fragmentation efficiency and relates the degree of unfolding and ECD efficiency to protein structure in the gas phase.

Novel Aspect

In-source ion activation coupled with IM-MS and ECD ion fragmentation for greatly improved sequence coverage and native protein structural analyses.

References

1. Eschweiler, J. D. ; Rabuck-Gibbons, J. N. ; Tian, Y. ; Ruotolo, B. T. Anal. Chem. 87, 11516 (2015).

ThP-169 / ANALYSIS OF INSULIN AND INSULIN ANALOGUES FROM DRIED BLOOD SPOTS BY MEANS OF LC-IM-HRMS

Andreas Thomas - Mario Thevis

German Sport University Cologne, Biochemistry, Cologne

Keywords: Ion mobility, high resolution MS, insulin, dried blood spots

Introduction

Analyzing peptides such as insulin from dried blood spots DBS represents an urgent task in various fields of analytical chemistry (e.g. doping controls, forensic- or clinical chemistry) due to simplified sampling, transfer and storage. Coupling liquid chromatography, ion mobility separation with high resolution mass spectrometry enables the sensitive detection of insulin and its synthetic analogues from DBS at physiologically relevant levels.

Methods

The blood spots (20 µL) were extracted from pretreated DMPK cards, immunopurified with magnetic beads and anti-insulin antibodies, separated by reversed-phase liquid chromatography (one-pump-trapping) and detected by high resolution mass spectrometry (time-of-flight) including an ion-mobility separation (T-wave IMS). The method was fully validated considering qualitative result interpretation.

Results

The evaluation of the diagnostic ion traces were performed under consideration of the respective retention times, accurate masses and specific drift times for each of the target peptides. With this strategy, validation of the method yielded specific and sensitive (< 0.5 ng/mL) results for human insulin and all relevant synthetic analogues. Recoveries range at approximately 30 % with precise (CV<25%) and linear results within the working range of the method (0.5 to 10 ng/mL). The proof-of-principle was shown by further analyzing DBS samples from healthy volunteers at non-fasting level (endogenous insulin detected only) and one diabetic volunteer (Type I). Here the rapid acting synthetic insulin Novolog was clearly detected. Generally, in this study was shown that additional separation by ion mobility combined with the accurate mass enhances the results significantly.

Conclusions

Analyzing insulins from urine or blood samples represents an established approach in doping controls [1].

The present study described now a fit-for-purpose method to analyze insulin from DBS cards. The application of DBS sampling becomes more and more interesting in various analytical fields. Here, the combination of immunoextraction, chromatography and high resolution MS with ion mobility separation offers the required sensitivity (sub ng/mL level).

Novel Aspect

To best of our knowledge the study shows the first assay to determine insulin and its analogues from dried blood spots in physiologically relevant levels with LC-MS.

References

1. Thomas A, Schänzer W, Thevis M., Drug Test Anal., 6, 1125-32 (2014)

ThP-170 / 5 HOUR SINGLE SHOT HUMAN PROTEOME ANALYSIS USING LC-FAIMS-MS/MS

Alexander Hebert⁽¹⁾ - *Satendra Prasad*⁽²⁾ - *Michael Belford*⁽²⁾ - *Derek J Bailey*⁽²⁾ - *Susan Abbatiello*⁽²⁾ - *Romain Huguet*⁽²⁾ - *Graeme McAlister*⁽²⁾ - *Elroy Wouters*⁽²⁾ - *Jean-Jacques Dunyach*⁽²⁾ - *Michael Westphall*⁽¹⁾ - *Joshua Coon*^(1,3)

University of Wisconsin-Madison, Biotechnology Center, Madison (1) - Thermo Fisher Scientific, Research and Development, San Jose (2) - University of Wisconsin-Madison, Department of Biomolecular Chemistry, Department of Chemistry, Morgridge Institute for Research, Madison (3)

Keywords: Proteomics, FAIMS, Ion Mobility, Quantitation, and FTMS

Introduction

Limitations in separation peak capacity and dynamic range restrict identifications in single shot proteomics experiments to ~ 5,000 protein groups. Pre-fractionation can skirt these issues; however, it comes at the cost of increased labor and number of samples for MS. Here we demonstrate that High Field Asymmetric Waveform Ion Mobility Spectrometry (FAIMS) coupled with an Orbitrap MS can allow near whole proteome sequencing without the need for pre-fractionation.

Methods

Proteins were extracted from a human cells (K562) and digested with trypsin. Peptides were loaded onto a nanoflow capillary column (30cm 1.7µm BEH C18), ionized, and further separated with a cylindrical electrode FAIMS device prior to injection into an Orbitrap Lumos. The FAIMS system operated with a 3MHZ bi-sinusoidal waveform (5000V peak to peak), electrode temperature of 100°C, and CVs between -35 and -110V dc. Protein identifications were filtered to 1% FDR.

Results

We describe a compact, cylindrical, helium free, rapid CV stepping, FAIMS source that easily interfaces between the ESI emitter and the MS inlet. This device induces an aerodynamic Coanda effect that entrains and diverts ions away from the inner electrode surface, yielding a high ion transmission efficiency (>50%). 20 ms transit times were obtained by minimizing the ion separation gap and electrode length (2.5mm to 1.5mm and 24mm to 15mm, respectively).

Performance test for single 90 minute analyses with and without FAIMS yielded 6,312 and 5,267 proteins, respectively. Triplicate FAIMS analysis of the same sample with 3 separate CV values (-40V, -60V, -80V) yielded 87,527 unique peptides and 7,427 proteins. To obtain similar performance without FAIMS offline high pH reversed phase chromatography (3 fractions total) had to be performed (95,745 unique peptides and 7,433 proteins). We conclude that FAIMS can accomplish gas-phase separations that offer comparable proteomic depth to conventional two-dimensional approaches.

Conclusions

Beyond these results we present a detailed characterization the z , m/z , mass and sample purity dependence of FAIMS. We compare quantitative (LFQ) reproducibility and identification overlap with and without FAIMS. Lastly we explore the possibility of analyzing with human proteome in a single-shot analysis where we identified 8,023 proteins with intra-analysis CV switching across 4 CVs (-45V, -60V, -75V, and -90V) during a single 5 hour experiment.

Novel Aspect

A novel FAIMS device enabling gas-phase fractionation as a viable alternative for rapid and deep proteomic analysis.

ThP-171 / DIFFERENTIATION OF PANAX GINSENG, PANAX QUINQUEFOLIUS, AND GINSENG PRODUCTS BY USING DIFFERENTIAL ION MOBILITY SPECTROMETRY COUPLED WITH TANDEM MASS SPECTROMETRY

Ri Wu - Tak Wah Dominic Chan

Faculty of Science, The Chinese University of Hong Kong, Department of Chemistry, N. T.

Keywords: Differential Ion Mobility Spectrometry, Tandem Mass Spectrometry, Ginseng, Collision-Induced Dissociation, Electron-Induced Dissociation

Introduction

Ginseng are highly treasured medicinal herb in traditional Chinese medicines. However, authentication of the origins of the ginseng products is not an easy task [1]. An DMS-MS/MS based method is established for the differentiation of Panax ginseng (Oriental ginseng) and Panax quinquefolius (American ginseng).

Methods

This method is based on the separation of ginsenosides present in the ginseng extracts using Differential Ion Mobility Spectrometry (DMS), followed by data dependent acquisition with Tandem Mass Spectrometry. Collision-induced dissociation (CID) and electron-induced dissociation (EID) were applied to generate diagnostic fragment ions for the differentiation of structural isomers [2, 3].

Results

The presence of ginsenoside Rf and pseudoginsenoside F11 in the Oriental and American ginsengs was selected as the criteria of differentiation. A good separation of isomeric ginsenosides was achieved in the Differential Ion Mobility Spectrometry using aliphatic alcohols as polar gas-phase modifiers. Data Dependent Acquisition results of DMS coupled with Collision-Induced Dissociation Mass Spectrometry was highly efficient compared with LC-MS/MS. Continuous ion separation in DMS chamber enabled the utilization of other dissociation methods that required longer acquisition time, such as EID. In comparison with collision-induced dissociation (CID), EID provides more structural informative ions and eliminates glycosyl group migration of Rg1 and Rf.

Conclusions

The spatial separation of isomeric ginsenosides in DMS and the generation of structural diagnostic ions for isomer differentiation by EID could eliminate ambiguities in differentiating ginseng and ginseng products. This study demonstrates a fast and effective method for differentiation of Oriental ginseng and American ginseng; and for characterization of commercial ginseng products.

Novel Aspect

High performance DMS-MS/MS method for differentiation and characterization of ginseng and ginseng products.

References

1. T. W. D. Chan, P. P. H. But, S. W. Cheng, I. M. Y. Kwok, F. W. Lau, H. X. Xu, *Anal. Chem.*, 72, 1281-1287 (2000).
2. Y. -L. Elaine W., Xiangfeng C., Wan L., Ze W., Y. -L. Winnie H., Ri W., T. -W. Dominic C., *Anal. Chem.*, 88, 5590-5594 (2016).
3. Xiangfeng C., Ze W., Y. -L. Elaine W., Ri W., Feng Z., T. -W. Dominic C., *Mass Spec. Rev.*, 1-18 (2018)

ThP-172 / ULTRA-HIGH THROUGHPUT SCREENING FOR ENGINEERED SYNTHETIC BIOLOGY LIBRARIES BASED ON SAMPLE CHIPS COMBINED WITH DESI MSIMAGING

Cunyu Yan⁽¹⁾ - **Eva Garcia-Ruiz**⁽¹⁾ - **Mark Dunston**⁽¹⁾ - **Katherine Hollywood**⁽¹⁾ - **Neil Swainston**⁽¹⁾ - **Adrian Jervis**⁽¹⁾ - **Emrys Jones**⁽²⁾ - **Aurore Lejeune-Dodge**⁽³⁾ - **Emily Kempa**⁽¹⁾ - **Yizhi Cai**⁽¹⁾ - **Eriko Takano**⁽¹⁾ - **Nigel Scrutton**⁽¹⁾ - **Sabine Flitsch**⁽¹⁾ - **Perdita Barran**⁽¹⁾

Manchester Institute of Biotechnology, University of Manchester, Manchester (1) - Waters Corp., Waters, Wilmslow (2) - Labcyte Europe, Labcyte Europe, Staffordshire (3)

Keywords: DESI, mass spectrometry, ion mobility, high throughput, synthetic biology

Introduction

The screening of variant libraries is an important yet costly and time consuming step in full Design-Build-Test-Learn pipelines[1] central to many current synthetic biology strategies. By employing liquid handling robotic platforms for automated sample-chip preparation, our new HTP strategy can generate visualized imaging data on the scale of 12 sample/min with only nL quantities needed for each sample.

Methods

MSImaging was performed by using a DESI 2D Omni Spray Ion Source combined with Synapt G2 Si HDMS mass spectrometer[2]. Sample chips were prepared using an Echo 550 liquid handling robotics system and cell strain supernatants transferred onto glass slides covered with a nylon membrane. The maximum capacity of each sample chip was 1536 samples/plate. The DESI imaging pixel size was set as 100µm X 100µm, with a scan rate of 15ms/spectrum. The acquired data were processed with HDImaging software.

Results

(i) Ambient and direct sampling conditions dramatically reduce liquid handling of samples; (ii) target compound titer measurements are label-free, and can thus be applied to different production pathways for the screening of intermediates and final products from synthetic biological libraries; (iii) the robotic sample chip preparation and low volumes required allow time course tracking and real time monitoring; (iv) the HTP sensitivity can reach µM level (*e.g.* in the case of pinocembrin flavonoid screening), which is far lower than the normal titres, and demonstrates the potential to apply this approach to other engineered pathways; (v) higher sensitivity is available when triple quadrupole mass spectrometry is applied; (vi) visualized MSImaging measurement provides highly accessible data for upstream Learning and subsequent re-Design-Build.

Conclusions

Our super-high throughput strategy, based on large capacity sample chips combined with MS imaging technology, makes medium and large scale directed evolution synthetic biology library screening possible over reasonable time scales. The superfast measurement data generated can provide sufficient data to allow the application of machine learning to model and debug malfunctioning synthetic biology pathways.

Novel Aspect

The super-high throughput MS screening method combines the automated production of sample chips and

MS imaging that presents a powerful new capability for the screening of large synthetic biology libraries.

References

1. Carbonell P., Jervis A., Robinson C., et al, Nat. Commun. Bio., 1, in press (2018).
2. Yan C., Parmeggiani F., Jones E., et al, J. Am. Chem. Soc., 139, 1408–1411 (2017)

ThP-173 / STRUCTURAL INSIGHTS INTO POTENTIALLY TOXIC CONFORMATIONS OF β -AMYLOID (1-42) OLIGOMERS

Rene Zangl - Janosch Martin - Tobias Lieblein - Nina Morgner

Goethe University, Institute of Physical and Theoretical Chemistry, Frankfurt

Keywords: Alzheimer's disease, Beta-Amyloid (1-42), Native-MS, LILBID-MS, ESI-IM-MS

Introduction

β -Amyloid (1-42) ($A\beta_{42}$), the protein involved in Alzheimer's disease (AD) has been demonstrated to generate β -sheet forming oligomers leading to fibrils and plaques. Whereas fibrils have long been suspected to cause neurodegenerative disorders, recent studies suggest misfolded oligomers as the origin of neurotoxic effects in AD. Hence the early oligomerization stage stands in focus of inhibitor development to suppress the formation of toxic oligomers. [1]

Methods

Additionally to laser induced liquid bead ion desorption mass spectrometry (LILBID-MS), which gives an insight in the early aggregation kinetics of $A\beta_{42}$ in absence and in presence of aggregation inhibitors, electrospray ionization mass spectrometry coupled with ion mobility (ESI-IM-MS) is used to reveal the inhibitor's structural influences on $A\beta_{42}$. [2]

Results

Whereas LILBID-MS allows to observe the aggregation kinetics and demonstrates the effect of well-known and novel aggregation inhibitors at the oligomeric level ESI-IM-MS provides a detailed look on several inhibitor binding events. Distinct m/z signals reveal the presence of several gas-phase conformations in their arrival time distribution (ATD) of the same oligomer order. Due to inhibitor interactions, the corresponding inhibitor bound states differ in their collision cross sections (Ω), as well in their gas-phase stabilities. Those indicate structural rearrangements which thereby affect distinct aggregation pathways differing in their kinetic behavior.

Conclusions

The investigation of inhibitors by ESI-IM-MS provides an insight into binding behavior of $A\beta_{42}$ at the oligomeric level. A closer look at inhibitor bound states reveals differences due to conformational changes which potentially lead to non-toxic oligomeric states. Estimating aggregation kinetics of $A\beta_{42}$ in absence and in presence of inhibitors by LILBID-MS additionally to the findings by ESI-IM-MS reveals the structural influence on the aggregation behavior.

Novel Aspect

The combination of LILBID and ESI-IM-MS provides promising results in understanding the aggregation of self-assembling peptides and their prevention by inhibitors.

References

1. Woods L. A., Radford S. E., Ashcroft A. E., Biochimica et Biophysica Acta (BBA) - Proteins and Proteomics, 1834, 1257-1268 (2013).

2. Stark T., Lieblein T., Pohland M., Kalden E., Freund P., Zangl R., Grewal R., Heilemann M., Eckert G. P., Morgner N., Göbel M. W., *Biochemistry*, 56, 4840-4849 (2017).

ThP-174 / TIMSTOF PRO AND PASEF MULTIPLYING SEQUENCING SPEED AND SENSITIVITY IN MASS SPECTROMETRY-BASED PROTEOMICS

Andreas-David Brunner⁽¹⁾ - *Florian Meier*⁽¹⁾ - *Scarlet Beck*⁽²⁾ - *Heiner Koch*⁽²⁾ - *Markus Lubeck*⁽²⁾ - *Niels Goedecke*⁽²⁾ - *Michael Krause*⁽²⁾ - *Nicole Drechsler*⁽²⁾ - *Jürgen Cox*⁽³⁾ - *Oliver Raether*⁽²⁾ - *Matthias Mann*⁽¹⁾

Max Planck Institute of Biochemistry, Department for Proteomics and Signal transduction, Planegg, Germany (1) - Bruker Daltonik GmbH, R&D, Bremen, Germany (2) - Max Planck Institute of Biochemistry, Systems Biochemistry, Planegg, Germany (3)

Keywords: Trapped ion mobility spectrometry, High sensitivity, High speed, Proteomics

Introduction

Applications of proteomics technologies to biomedical research require further developments of mass spectrometry (MS) technology to overcome longstanding limitations in speed, sensitivity, and robustness. Trapped ion mobility spectrometry (TIMS) coupled to a quadrupole time-of-flight (qTOF) analyzer in combination with the 'Parallel Accumulation – SErial Fragmentation' (PASEF) scan mode has shown much promise in this regard. [1]

Methods

TIMS separates ions based on their mobility and releases them sequentially into the MS. PASEF synchronizes the MS2 precursor selection with the TIMS elution time - multiplying sequencing speed and sensitivity by multiple precursor selection per TIMS scan [2]. An algorithm that detects precursors in 4 dimensions was developed. Data analysis is implemented in MaxQuant. We evaluated the TIMS-TOF Pro performance with regard to sensitivity, speed, and reproducibility.

Results

In 2h single-shot measurements of 200ng HeLa digest, our setup detected more than 350,000 multiply charged features and identified more than 6,000 protein groups with an high quantitative reproducibility between replicates ($R2 = 0.97$) and a scan speed of >120 Hz. From only 10 ng of HeLa digest, which converts to approximately 30 HeLa cells, we identified more than 2.400 protein groups in 60 min and more than 2.450 protein groups in only 30 min, highlighting the sensitivity and sequencing speed of our setup. [3] We demonstrate the combination of the exceptional scan speed and high sensitivity of our mass spectrometry setup for turbo protein-protein interactome studies in yeast pulldown experiments. We identify all known interaction partners from protein complexes distributed across the whole protein copy number range and from only 50-100 ng of input material in 23 min measurements.

Conclusions

The additional separation of co eluting ions and signal compression in sharp ion mobility peaks in the additional TIMS dimension increases sensitivity and dynamic range. This allows high speed and high sensitivity applications, and routine large-scale measurements of peptide collisional cross sections (CCS), reflecting their 3D structure, which adds a new and yet underexplored level of information to conventional bottom-up proteomics workflows.

Novel Aspect

Robust, high sensitivity and high speed shotgun proteomics on a novel TIMS-qTOF platform combined with the PASEF scan mode.

References

1. Ridgeway, ..., Park, Int. J. Mass Spec., 425 (2018)
2. Meier, ..., Mann, J. Proteome Res., 14, 5378-5387 (2015)
3. Meier, ..., Mann, 2018, bioRxiv (2018)

ThP-175 / CROSS SECTION OF SELECTED EXTRACTABLE AND LEACHABLE COMPOUNDS: INFLUENCE OF MOBILE PHASE COMPOSITION

Vincent Hanot⁽¹⁾ - Rodrigo Feliciano⁽¹⁾ - Peifeng Hu⁽²⁾ - Christopher Jones⁽²⁾

Baxter, Mass Spectrometry, Braine L'alleud, Belgium (1) - Baxter, Mass Spectrometry, Round Lake, United States (2)

Keywords: Ion Mobility, CCS, E&L, Identification, Mobile Phase

Introduction

Ion mobility coupled to mass spectrometry (IM-MS) allows the separation of molecular species based on their mass, charge and shape, providing an additional separation dimension to chromatography. In this work, we compare the reproducibility of collision-cross section (CCS) measurements for analytes separated using different liquid chromatography (LC) mobile phase conditions to evaluate the robustness of LC-IM-MS for extractables identification.

Methods

In order to study the impact of pH on adduct formation of known impurities from pharmaceutical packaging material by IM-MS, several aqueous mobile phases with different pH were evaluated. IM-MS detection was performed on a Waters VION in positive and negative modes between m/z 80-1000. Each mixture of standards analysed in triplicate was bracketed by QC analysis ($n=5$) and CCS values for each compound were collected.

Results

The reproducibility of CCS values has been measured for analytes originating from different LC mobile phases. Different pH were evaluated (10 mM ammonium acetate + 0.05% CH₃COOH, 0.1% HCOOH and 0.001% HCOOH) using methanol as organic mobile phase during a 35 min gradient method. 10 μ L of standards or QC samples were injected into an Agilent Zorbax Eclipse Plus, 2.1 x 100 mm at 40°C with a 0.5 mL/min flow. All molecules were detected in negative mode as [M-H]⁻. In positive mode, each analyte was detected as [M+H]⁺, [M+Na]⁺, [M+K]⁺, [M+NH₄]⁺, or different combinations of these adducts. The precision of CCS values was $\leq 1.13\%$ and $\leq 0.57\%$ for QC and samples, respectively. The difference between expected and observed CCS (Δ CCS) ranged between [-1.60; 1.02%] and [-1.94; 1.30%] for QC and samples, respectively. A very good precision and "accuracy" were measured independently of the pH of the mobile phase. The results of this study also showed that the relationship between the adduct formed and the measured CCS value is not straightforward.

Conclusions

These results showed that CCS values for small organic ions are pH-independent, unlike the retention time. The very good reproducibility of CCS values highlighted the usefulness of CCS values in identifying compounds and can serve as a complementary fingerprint for the identification of unknown compounds in E&L studies. In the identification process, the adduct formation should be considered as it will impact the measured CCS value.

Novel Aspect

This work highlights the relevance of using CCS values as an additional dimension to identify unknown compounds in E&L studies.

ThP-176 / LIPID AND FATTY ACID ISOMER ANALYSIS USING A TRAVELLING WAVE CYCLIC ION MOBILITY SEPARATOR

Mike McCullagh⁽¹⁾ - *Martin Palmer*⁽¹⁾ - *Giorgis Isaac*⁽²⁾ - *Russell Mortishire-Smith*⁽¹⁾ - *James Langridge*⁽¹⁾ - *Johannes Vissers*⁽¹⁾

Waters Corporation, Wilmslow, United Kingdom (1) - Waters Corporation, Milford, United States (2)

Keywords: IM, lipid, isomer, resolution, traveling-wave

Introduction

Lipids represent a diverse group of biomolecules that have an essential role in structural, storage and signaling processes in living systems. Class separation is readily achieved using chromatographic and MS based identification techniques; however, the analysis remains challenging due to the chemical structure diversity and isobaric nature of lipids. The addition of ion mobility (IM) to lipidomic workflows enhances separation and improves isomer resolution.

Methods

Data were collected on a hybrid quadrupole cyclic IM orthogonal acceleration time-of-flight (Q-cIM-oaToF) instrument. Ion mobility separation is achieved using a multi-pass travelling-wave cIM separator, where increasing the number of passes around the device allows of the increase in both mobility resolution and ion residence time. MS and CID fragmentation data were obtained on precursor IM separated lipids followed by ToF mass measurement.

Results

Unsaturated free fatty acid (FA) standards, differing in chain length and number of cis/trans configurations, were chosen to determine the degree of IM separation required to separate lipid isomers. FAs represent the simplest class of lipid components, yet they exist in the structure of each lipid categories. In all direct infusion cIM-MS measurements, FAs with cis-double bond orientations, introduced as two component mixtures, were found to be more compact than those with trans-orientations. Moreover, the cis- and trans-orientations for the monounsaturated FAs were distinguishable. A different number of cycles through the cIM separator, thereby increasing the effective path length/resolution, were required to achieve a similar degree of IM separation for mono-unsaturated FAs of differing chain length. Unsaturated FAs with two or more double bonds, separated by two mid-chain carbons, could not be distinguished. The separation of other lipid classes by cIM is currently under study.

Conclusions

A Q-cIM-oaToF geometry has been characterized and successfully applied to the IM separation of isomeric, mono-unsaturated FAs. The required IM resolution was found to be a function of chain length. Shorter, structurally more rigid and compact FAs were discriminated at reduced resolution, as could longer chain mono-unsaturated FAs as a result of partial chain back-folding. Following IM separation, isomeric lipids were successfully CID fragmented and identified.

Novel Aspect

High resolution ion mobility separation of cis/trans lipid isomers using a using a multi-pass travelling-wave cIM separator.

References

1. Giles K., Wildgoose J. L., Pringle S., Garside J., Carney P., Nixon P., Langridge D. J. In 62nd Annual ASMS Conference on Mass Spectrometry and Allied Topics, Baltimore, MD, June 15–19, 2014.

ThPS-S08 - Petroleomics, hydrocarbons and biofuels

Chairs: Alan G. Marshall, Giorgio Mellerio

ThP-177 / DETAILED LOOK IN SARA FRACTIONATION: SATURATES, NOT ONLY SATURATED COMPOUNDS!

Zahra Farmani - Wolfgang Schrader

Mass Spectrometry dept., Max-Planck-Institut für Kohlenforschung, Mülheim a. d. Ruhr

Keywords: SARA, Saturates, HRMS, direct analysis

Introduction

Saturates, Aromatic, Resins and Asphaltenes (SARA) fractionation is a decade old fractionation method of crude oils. Saturates, referring to their chemically descriptive name, are known to contain only aliphatic compounds. To date, there are very few saturates centered studies of which all were performed using GC techniques. Therefore, their direct MS analysis were neglected and never studied in detail.

Methods

SARA fractionation was performed on different crude oils using a modified IP-143 method. [1] Saturates were dissolved in Toluene: Methanol (1:1) to the final concentration of 250 ppm. results of the direct injection HRMS analysis by ESI-(+)-FTMS for all of the samples were compared. In order to get further comparable information, all of the samples were further analyzed using HRGC-EI-MS on a GC-Orbitrap MS.

Results

Direct injection HRMS was successfully performed for all of the samples. Comparison of direct injection analysis shows that saturates of different oils, poses different chemical characteristics including their heteroatom containing classes and the length of their aliphatic chains. results of the HRMS analyses are contradicting the common knowledge about the characteristics of saturates. In this regards, the compounds found in this fraction do not necessarily only contain aliphatic chains and saturated rings, but also many olefinic chains and unsaturated rings. More than 15000 different elemental compositions could be calculated from the data, showing a much wider variety of compounds than expected from a saturate fraction. The GC-MS measurement results show presence of volatile and aliphatic compounds that are obviously only a part of the saturates fraction, which are detectable using GC technique.

Conclusions

Results of the HRMS analysis prove that the characteristics of saturates fraction differ from the contemporary knowledge consisting of only aliphatic compounds. This is because SARA fractionation is solubility, polarity and adsorption based separation method which cannot give a clear separation and results overlaps between fractions.

Novel Aspect

Saturates fraction contains not only GC detectable aliphatic compounds, but many olefinic and aromatic compounds which can be analyzed by direct injection HRMS.

References

1. H. Alboudwarej, J. Beck, W. Y. Svrcek, H. W. Yarranton, K. Akbarzadeh., Energy Fuels, 16, 462-469 (2002).

ThP-178 / EVALUATION OF DEGRADATION DEGREE OF LUBRICATING OIL BY THERMAL DESORPTION AND PYROLYSIS COMBINED WITH DART-MS (TDP/DART-MS)

Chikako Takei - Yukiko Otake - Kenichi Yoshizawa

BioChromato, Inc., Application support, Fujisawa, Japan

Keywords: Lubricating oils, evaluation of degradation degree, TDP/DART-MS

Introduction

Lubricating oils are composed of a base oil and various additives. Generally, since complicated pre-treatment is necessary prior to analysis, this extending the time and increases the effort required to obtain the final results. Thus, analysis without pre-treatment would be desirable for research and development. Therefore, the purpose of this work is to directly and rapidly evaluate the degradation degree of lubricating oil.

Methods

Automobile lubricating oil, virgin and after 1000 or 5000 km driving were used. Thermal desorption and pyrolysis (TDP) combined with DART-MS was composed of Q-TOF mass spec equipped with DART ion source and ionRocket TDP device. The samples were used to ~1 μ L and put into the sample pot. TDP temperature condition was room temperature to 600°C, a rate of 100°C/min. To structurally and compositionally characterize samples, Kendrick Mass Defect (KMD) analysis was used.

Results

The mass spectra were detected for the thermal desorption and pyrolysis products of all samples. In the mass spectra detected at 300°C, additives such as a phenol-type antioxidant, an amine-type antioxidant, and a salicylic acid-type detergent dispersant were detected in all the samples. There was no significant difference in the ratio of additive compounds at 0 and after 1000 km of driving, but after 5000 km driving, the concentrations of the salicylic acid-type detergent and phenolic antioxidant were remarkably decreased. In the mass spectra recorded at 450°C the base oil dominated the spectrum at 0 km, while after 1000 km of driving, glycol compounds were also detected, which increased in concentration after 5000 km of driving. Additionally, using KMD plots it differences among the samples were clearly observed. By setting the Kendrick Mass unit based on the mass of the repeating structure, compounds with common repeat units aligned in the horizontal direction on the KMD plot, whereas compounds with different structures shifted.

Conclusions

TDP/DART-MS enables the detection of additives and base oil in lubricating oils by gradient heating without pre-treatment and can be used to evaluate the degradation extent of lubricating oils quickly and easily by monitoring the amount of antioxidants as a degradation marker. Moreover, a combination of TDP/DART-MS and KMD analyses enabled clear recognition of slightly differences in the base oils between the virgin and degraded samples. Thus, the method described herein can be used for the evaluation of the degree of degradation of lubricating oils.

Novel Aspect

A combination of TDP/DART-MS and KMD analysis is useful for evaluation of degradation degree of lubricating oils.

ThP-179 / RADICALS GENERATED BY NANOPARTICLES IN CONTACT WITH LIQUID WATER USING SPIN TRAPPING AND MASS SPECTROMETRY

Per Axel Clausen⁽¹⁾ - Vivi Kofoed-Sørensen⁽¹⁾ - Asger W. Nørgaard⁽²⁾ - Nicklas R. Jacobsen⁽¹⁾ - Keld A. Jensen⁽¹⁾

National Research Centre for the Working Environment, New Technologies, Copenhagen Ø (1) - Novo Nordisk, Analytical Department, Bagsvaerd (2)

Keywords: Nanoparticles, radical generation, spin trapping, MS

Introduction

Inhalation of particles may cause airway inflammation induced by reactive oxygen species (ROS) generated by the particles. ROS are often measured by non-specific methods *e.g.* using 2',7'-dichlorodihydrofluorescein (DCFH₂) and comprises both radicals (hydroxyl, superoxide) and non-radicals (hydrogen peroxide) which both contribute to the response of DCFH₂. The aim of the study [1] was to measure specific radicals formed by nanoparticles in contact with water.

Methods

The method is based on spin trapping using diethoxyphosphoryl-5-methyl-1-pyrroline (DEPMPO) that traps radicals right after their formation. In short, nanoparticles (NPs) were dispersed in pure water with added spin trap, incubated over-night, and centrifuged. The supernatant was analysed with LC-MS and the "pellet" consisting of the NPs was analysed with MALDI-TOF-MS to directly measure the spin trap radical products adsorbed to the surface of the NPs.

Results

It turned out that the major part of the spin trap radical products was adsorbed to the particle surfaces. Carbon black (Printex 90) showed by far the strongest radical generation as compared to silica NP (NM-200), with carbon nanotubes (NM-400) in between. The preliminary results showed that the Printex 90 produced mainly the hydroxyl radical ($\cdot\text{OH}$) but also the hydroperoxyl radical ($\cdot\text{OOH}$). NM-400 produced only $\cdot\text{OH}$ and NM-203 produced no radicals. The total radical formation potential (sum of $\cdot\text{OH}$ and $\cdot\text{OOH}$ from both supernatant and "pellet") of Printex 90 was 13-46 nmol/mg NP and 8-20 nmol/mg NP of NM-400 depending on the amount of NP dispersed in the water.

Conclusions

When studying the specific radical formation of NPs using spin traps and MS techniques it is necessary to be aware, that a substantial part of the spin trap radical products may be adsorbed to the surface of the NPs and has to be measured independently. The preliminary results show that radical formation by different types of NP may range from no radical formation of silica NP and up to about 50 nmol/mg NP of carbon black. All mechanisms behind ROS production are not known.

Novel Aspect

Quantitation of radicals formed by nanoparticles in their dispersion media and on their surface using spin trapping and MS has to our knowledge not been reported previously.

References

1. The present study was part of the Danish NanoSafety Centre.

ThP-180 / A NOVEL NANOCOMPOSITE COATING BASED ON SOL-GEL TITANIA/HYDROXYAPATITE FOR SOLID-PHASE MICROEXTRACTION COUPLED TO GC-MS

Fereidoon Rahmani⁽¹⁾ - *Ali Es-haghi*⁽²⁾ - *Mohammad-Reza Milani Hosseini*⁽¹⁾ - *Afsaneh Mollahosseini*⁽¹⁾

Iran University of Science and Technology, Department of Analytical Chemistry, Tehran (1) - Razi Vaccine & Serum Research Institute, Agricultural Research, Education and Extension Organization (AREEO), Department of Physico Chemistry, Karaj (2)

Keywords: Solid-phase microextraction; Sol-gel technology; Hydroxyapatite; Titania; GC-MS.

Introduction

The underlying principle of SPME technique is based on the analyte distribution between the sample matrix and the extraction phase [1]. Titania and hydroxyapatite (HAP) have received great interest in analytical chemistry due to very good physical and chemical properties [2]. HAP can improve the surface area, grain size and thermal stability of TiO₂. Addition of TiO₂ changes the surface morphology as well as an increase of the mechanical and chemical stability of HAP [3].

Methods

A non-aqueous organic/inorganic sol-gel rout was applied for synthesis of HAP sol and then a transparent stable titania sol was prepared. After repeating steps including; dipping, drying and sintering, a sol-gel coating was formed on the modified stainless steel (SS) wire. Finally, applicability of the prepared fiber was assessed for the HS-SPME of BTEX as model analytes from water samples followed by GC-MS.

Results

The thickness, morphology and composition of the prepared sol-gel coatings; TiO₂, HAP and HAP/TiO₂, was investigated and analyzed by using of FESEM and EDX. The coating thickness was estimated below 1μm when the coating process was performed 5 times. This is much less than the thickness of the thinnest conventional coating. The FESEM images showed a fully cracked and quite non-uniform coatings in terms of porosity and particle size by use of pure TiO₂ or HAP. Addition of HAP sol to the pure TiO₂ sol, considerably changes the morphology of HAP/TiO₂ coating; the coating is perfectly uniform and the size of the particles in HAP/TiO₂ nanostructure is increased compared to the pure TiO₂ coating. Nanoparticles (<50 nm) and pores (less than 100 nm) have been spread uniformly throughout the wire, providing a high surface area and a high sample loading capacity.

The obtained chromatograms of BTEX after HS-SPME by prepared TiO₂, HAP and HAP/TiO₂ coatings in optimum extraction conditions testify that HAP/TiO₂ coating exhibited high extraction capacity and have better peak shapes.

Conclusions

The proposed sol-gel method for the synthesis of HAP/TiO₂ and the coating process on SS wire are very cost-effective, simple, fast and structurally controllable. The obtained chromatograms of BTEX after HS-SPME showed that HAP/TiO₂ coating exhibited high extraction capacity. Very LOD and LOQ values were obtained for BETX.

Novel Aspect

A novel titania/hydroxy apatite nanocomposite coating fiber was prepared for first time through a relatively fast sol-gel technique.

References

1. Pawliszyn J., Handbook of solid phase microextraction, 1st edition, Elsevier Inc, Waltham, MA, USA. (2012).
2. Kim T. Y., Alhooshani Kh., Kabir A., Fries D. P., Malik A., J. Chromatogr. A, 1047, 165–174 (2004).
3. Gshalaev, V. S., Demirchan, A. C., Hydroxyapatite synthesis, properties and applications, Nova Science

Publishers, Inc., New York (2012).

ThP-181 / MASS ANALYSIS OF SURFACE IMMOBILIZATIONS OF PHENOL-CONTAINING MOLECULES THROUGH CHEMICAL MODIFICATIONS

Hye-Jeong Hwang⁽¹⁾ - *Inseoung Choi*⁽¹⁾ - *Young-Kwan Kim*⁽²⁾ - *Woon-Seok Yeo*⁽¹⁾

Konkuk University, Department of Bioscience and Biotechnology, Seoul (1) - Korea Institute of Science and Technology, Institute of Advanced Composite Materials, Carbon Composite Materials Research Center, Wanju-gun (2)

Keywords: Immobilization, Mannich reaction, Phenol-containing molecules, Self-assembled monolayers, Surface chemistry

Introduction

To construct applicable materials incorporated with biologically active phenol-containing molecules, strategies for immobilization of phenolic molecules on materials are required [1]. Although a number of immobilization methods have been accomplished and reported, mostly harnessing the phenol functionality, the practical use as a general immobilization method has been hampered due to their complicated chemical reactions and low reaction yield on surfaces [2-3].

Methods

Chemical functional groups, carboxylic acid, thiol, and azide, were introduced to phenolic molecules using Mannich reaction followed by immobilization to self-assembled monolayers on gold via various surface chemistries, carbodiimide coupling reaction, Michael addition, and 'click' reaction. As a practical application, we performed sulfotransferase enzyme assay on immobilized 2-naphthol [4]. The surface reactions were analyzed by MALDI-TOF MS and XPS.

Results

Naphthol, tyrosine, and hydroxyflavone were modified by Mannich reaction to possess carboxylic acid, thiol, and azide group, which were confirmed by MALDI-TOF MS. These functionalized phenolic molecules were then immobilized to the monolayers on gold that presented amine, maleimide, and alkyne via carbodiimide coupling reaction, Michael addition, and 'click' reaction, respectively. The immobilization of the synthesized molecules to the monolayers were analyzed by MALDI-TOF MS and XPS. The mass spectra clearly showed that all the surface conjugation reactions proceeded as expected, although the fragmented final products were observed. In addition, comparison of C1sXPS indicated that the molecules were immobilized on the SAMs successfully. As a practical application, the 2-naphthol immobilized monolayers were exposed to sulfotransferase SULT1A1 and 3'-phosphoadenosine 5'-phosphosulfate, which yielded 2-naphthol sulfate.

Conclusions

we describe a simple, fast, and reliable method for the surface immobilization of phenol-containing molecules such as flavonoids by introducing chemical functional groups using Mannich reaction while maintaining the hydroxyl functionality so that their biological activities were not compromised.

Novel Aspect

Our method can be a general and practical platform for the immobilization method of various phenol-containing molecules on surfaces of various materials.

References

1. Balasundram N., Sundram K., Samman S., *Oxid. Med. Cell Longev.*, 5, 270–278 (2009).
2. Arango D., Morohashi K., Yilmaz A., Kuramochi K., Parihar A., Brahimaj B., Grotewold E., Doseff A.I.,

- PNAS, 110, E2153-E2162 (2013).
- Lu Y., Irani N., Grotewold E., BMC Chemical Biology, 5, 3-12 (2005).
 - K. Halake, M. Birajdar, J. Lee, Journal of Industrial and Engineering Chemistry 35 (2016) 1–7
 - Falany C., Ström P., Swedmark S., Br. J. Clin. Pharmacol., 60, 632-640 (2005).

ThP-182 / NEW BINDING PROPERTIES OF PYRIDINE[4]ARENE CAPSULES STUDIED BY IM-MS AND IRMPD

Anniina Kiesilä⁽¹⁾ - *Anneli Kruve-Viil*⁽²⁾ - *Christoph A. Schalley*⁽²⁾ - *Perdita Barran*⁽³⁾ - *Jani Moilanen*⁽¹⁾ - *Elina Kalenius*⁽¹⁾

University Of Jyväskylä, Department of Chemistry, Jyväskylä (1) - Freie Universität Berlin, Institut für Chemie und Biochemie, Berlin (2) - University of Manchester, Manchester Institute of Biotechnology, Manchester (3)

Keywords: Supramolecular chemistry, Pyridine[4]arenes, Ion mobility mass spectrometry, IRMPD

Introduction

Pyridine[4]arenes [1] are macrocycles, which form dimeric and hexameric capsules via hydrogen bonding [2]. This type of capsules are widely studied in supramolecular chemistry as hosts, which can bind guest molecules inside the capsular cavity. We have recently discovered new binding properties of pyridinearene capsule, which was previously assumed to encapsulate only anionic guests [3-5]. The structures of new complexes have been studied using IRMPD and IM-MS.

Methods

Complexation of pyridine[4]arene with cationic guests was studied by ESI-Q-TOF MS (QSTAR from ABSciex). The structural features of observed complexes were studied in IRMPD experiments (performed with Varian FT-ICR) and DTIM-MS. DFT calculations have been used to visualize the gas phase structures.

Results

Complexation properties of pyridinearene were tested with different alkyl ammonium guests. In (+)ESIMS experiments, complexation was observed only with tetramethyl ammonium cation (TMA). In (-)ESI-MS spectra ternary complexes with TMA cation and two anions were observed to bound to pyridinearene dimer. These ternary complexes were studied in infrared multiphoton dissociation (IRMPD) and drift tube ion mobility experiments. Gained results supported the encapsulation of cationic guest while two anions were bound outside the capsule. The gas phase structures were visualized by DFT calculations.

Conclusions

New type of supramolecular ternary complexes have been observed in gas phase. Different mass spectrometric techniques have been utilized to study their structure. results were consistent and supported the encapsulation of cationic guest, while anions were bound outside the capsular cavity.

Novel Aspect

New type of supramolecular ternary complexes observed in gas phase. The structural features of these complexes were studied using IRMPD experiments and IM-MS.

References

- Gerkenmeier T., Mattay J., Näther C., Chem. Eur. J., 7 (2001).
- Evan-Salem T., Cohen Y., Chem. Eur. J., 13 (2007).
- Letzel M. C., Decker B., Rozhenko A. B., Schoeller W. W., Mattay J., J. Am. Chem. Soc., 126 (2004).
- Rozhenko A. B., Schoeller W. W., Letzel M. C., Decker B., Mattay J., New J. Chem., 37 (2013).
- Kiesilä A., Kivijärvi L., Beyeh K., Moilanen J., Groessl M., Rothe T., Götz S., Topić F., Rissanen K., Lützen A., Kalenius E., Angew. Chem. Int. Ed. 56, (2017)

ThP-183 / MASS SPECTROMETRIC ANALYSIS OF VICINAL DIOL-CONTAINING FLAVONOID MOLECULES VIA ENRICHMENT WITH BORONIC ACID-FUNCTIONALIZED PARTICLES

Eunjin Kim⁽¹⁾ - *Jihyeon Song*⁽²⁾ - *Hyejung Mok*⁽²⁾ - *Woon-Seok Yeo*⁽¹⁾

Department of Bioscience and Biotechnology, Bio/Molecular Informatics Center, Konkuk University, Seoul 05029 (Korea) (1) - Department of Bioscience and Biotechnology, Konkuk University, Seoul 05029 (Korea) (2)

Keywords: Boronic acid-functionalized particles (BAPs), Enrichment, Flavonoid, Matrix-assisted laser desorption/ionization time of flight mass spectrometry (MALDI-TOF MS)

Introduction

Detection and quantification of flavonoids which exhibit broad biological and pharmacological activity with little toxicity [1] has been hampered by their bad water-solubility and rapid metabolic elimination from human body [2]. This study concerns the development of an efficient enrichment method for facilitating the analysis of vicinal diol-containing flavonoid molecules using boronic acid-functionalized particles (BAPs) and MALDI-TOF MS.

Methods

Boronic acids on BAPs bind to vicinal diols to form boronate monoesters at basic pH. This complex remains intact during the enrichment process, and the vicinal diol-containing flavonoids are easily separated by centrifugation and subsequent acidic treatment. The selectivity and the limit of detection of our strategy were confirmed by MS analysis, and the validity was assessed by performing the detection and quantification of a model molecule quercetin in mouse organs.

Results

The selectivity of BAPs towards vicinal diol-containing flavonoids was evaluated in the presence of other flavonoids which have multiple hydroxyl groups and carbonyl groups that potentially interact with boronic acids. Next, we observed excellent enrichment feasibility for vicinal diol-containing flavonoids of our method in the presence of other flavonoids in highly excess amount. A calibration curve for quantification was constructed with various amounts of quercetin in the presence of an internal standard, and the limit of detection for quercetin was 75 pmol. The practical applicability was confirmed by quantifying quercetin spiked in fetal bovine serum with excellent accuracy and precision. Furthermore, we showed that our method can effectively enrich vicinal diol-containing flavonoids and facilitate the quantification of the targets present in animal tissues in low abundance. Finally, excellent recyclability of the BAPs with good reproducibility was confirmed indicating the practical usability and cost saving of our method.

Conclusions

We validated the selectivity and efficiency of BAPs towards vicinal diol-containing flavonoid molecules and constructed a calibration curve for quantifying quercetin. The feasibility and applicability of our method were assessed via the quantification of quercetin spiked in FBS and animal tissues. Furthermore, the recyclability of BAPs was confirmed. Taken together, this method is promising from the viewpoint of cost saving and practical applications.

Novel Aspect

Boronic acid functionality was harnessed for the analysis of flavonoids in spite of its low water-solubility and interference of other vicinal diol-containing biomolecules[3].

References

1. Tapas A. R., Sakarkar D., Kakde R., Trop. J. Pharm. Res., 7, 1089-1099 (2008).
2. Kim M. K., Park K. -s., Yeo W. -s., Choo H., Chong Y., Bioorg. Med. Chem., 17, 1164-1171 (2009).

3. Zhang L., Wang C., Wei Y., Biomed. Chromatogr., 29, 1133-1136 (2015).

ThP-184 / SPECTRAL INTERFERENCES AFFECTING THE NANOPARTICLES DETERMINATION USING SP-ICP-MS

Martin Loula - Antonin Kana - Oto Mestek

University of Chemistry and Technology Prague, Department of Analytical Chemistry, Prague

Keywords: single particle ICP-MS; arsenic nanoparticles; spectral interferences

Introduction

Analysis of nanoparticles (NPs) using single particle inductively coupled plasma mass spectrometry (sp-ICP-MS) is able to provide information on the size of NPs and their number concentration. However, the determination of some elements via ICP-MS is strongly hampered by spectral overlap from polyatomic species. Overcoming interferences of polyatomic ions, which can compromise results, is one of the challenges with ICP-MS.

Methods

The use of ICP-MS with reaction cell enables interference-free conditions to be obtained, even in the most demanding applications. This approach has also been proven to allow analysis of NPs in matrices containing interference precursors [1]. This work is focused on the determination of arsenic NPs in biological matrices because the arsenic suffers a polyatomic interference from ArCl^+ .

Results

The As NPs were prepared by borohydride-reduction of sodium arsenite according to previously published paper [2]. The use of NH_3 in reaction cell mode ($\text{FNH}_3=0.3$ mL/min) was demonstrated to be the most suitable to overcome the interference affecting trace As determination up to chloride concentration of 5.46 g/L, which corresponds to physiological saline solution. By using reaction cell, As NPs >50 nm can be detected. Because the sensitivity of signal was reduced not only by the reaction gas but also by salt content in sample, it was possible to provide accurate results only for As NPs >70 nm in at least 10 times diluted saline samples. In the case of <70 nm NPs, the histograms show incomplete distributions as result of the presence of the background signal.

Conclusions

sp-ICP-MS method utilizing reaction cell that allows determination of As NPs in presence of high chloride concentrations was introduced. Despite the reduction of sensitivity and therefore increased particle size detection limit, this method promises possibility of determining As NPs in biological matrices (*e.g.* blood serum).

Novel Aspect

sp-ICP-MS offers fast determination and characterization of As NPs. The applicability to biological matrices may have interesting uses in medicine.

References

1. Tharaud M., Gondikas A. P., Benedetti M. F., von der Kammer F., Hofmann T., Cornelis G., Journal of Analytical Atomic Spectroscopy, 32, 1400-1411 (2017).
2. Pal A., Saha S., Maji S. K., Kundu M., Kundu A., Advanced Material Letters, 3, 177-180 (2012).

Acknowledgements

Authors thank for financial support of the GACR project No. 17-00291S, specific university research (MSMT No 21-SVV/2018) and the project of institutional support of UCT Prague.

Shokoofeh Maghari - Parviz Ghezellou - Alireza Ghassempour

Shahid Beheshti University, Medicinal Plants and Drugs Research Institute Department: Phytochemistry, Tehran

Keywords: MALDI-TOF, mass spectrometry, nanomaterials, functionalized silica nanoparticles

Introduction

Nanomaterials have been widely used in laser desorption/ionization time-of-flight mass spectrometry as an alternative for matrix (LDI-MS) and as well for improving signal intensities in MALDI-TOF-MS [1]. In this regard, bare silica nanoparticles (SiO₂NPs) have been successfully utilized for the signal enhancement of peptides [2]. Here, we report influence of the two types of functionalized silica NPs in the development of MALDI-MS signals intensity of peptides.

Methods

The sulfobetaine silica NPs (SB-SiO₂-NPs) and amino propane sulfonate silica NPs (APS-SiO₂-NPs) were prepared according to the published procedure [3]. MALDI experiments were done by using the AB SCIEX4800 time-of-flight (TOF/TOF) mass spectrometer in reflectron positive mode, which is equipped with a 355 nm Nd:YAG laser. Typically, to acquire a single spectrum 800 laser shots were fired at a frequency of 200 Hz at 25% of the full laser power for samples.

Results

Among the different methods for sample preparation, two methods were selected for MALDI-TOF-MS experiment. In this regard, 1 μL of functionalized SiO₂-NPs was deposited onto the MALDI plate. Following that 1 μL mixture of sample solution, including peptide and matrix, was spotted on the top of dried silica NPs. In a parallel way we mixed 1 μL of functionalized SiO₂-NPs with 1 μL mixture of sample solution and after 3 min vortex, 1 μL of NPs/sample solution was directly spotted on to the MALDI plate and air dried for MALDI-MS analysis. In all analyses, the concentration of peptides in MALDI samples was maintaining constant. The maximum signal intensity for peptides was obtained when we used NPs as additive to the sample solution. SB-SiO₂-NPs had the most influence on analyte signal enhancement. Alteration the ratio of matrix and NPs in sample solution had a positive effect on signal enhancement. As well, concentrating sample with the coffee-ring effect during sample preparation, improved signals intensity.

Conclusions

This work demonstrates that functionalized SiO₂ NPs can be used for enhancement of ion signals in MALDI-TOF-MS analysis of peptides. This improvement on signals intensity is possibly related to the electrostatic interaction between NPs and peptides. Furthermore, alteration NPs surface charge have influence on this electrostatic interaction and change the amount of peptides ionization.

Novel Aspect

New functionalized silica nanoparticles have been utilized for the first time, since improving signals intensity of peptides in MALDI-MS analysis.

References

1. Yagnik G. B., Hansen R. L., Korte A. R., Reichert M. D., Vela J., Lee Y. J., *Anal. Chem.*, 88, 8926-8930 (2016).
2. Agrawal K., Wu H. F., *Rapid Commun. Mass Spectrom.*, 22, 283-290 (2008).
3. Knowles B. R., Wagner P., Maclaughlin S., Higgins M. J., Molino P. J., *ACS Appl. Mater. Interfaces*, 9, 18584-18594 (2017).

Martin B. Minameyer⁽¹⁾ - *Daniela Dasler*⁽²⁾ - *Ricarda A. Schäfer*⁽²⁾ - *Jakob F. Hitzenberger*⁽¹⁾ - *Frank Hauke*⁽²⁾ - *Thomas Drewello*⁽¹⁾ - *Andreas Hirsch*⁽²⁾

Friedrich-Alexander University (FAU), Physical Chemistry I, Erlangen-Nürnberg (1) - Friedrich-Alexander University (FAU), Organic Chemistry II, Erlangen-Nürnberg (2)

Keywords: LDI-ToF MS, Graphene, Porphyrin, Covalent Functionalization, Nanohybrid

Introduction

To characterize covalently functionalized graphene, several analytical methods are commonly applied. [1] Raman spectroscopy, for example, allows to detect the conversion of sp^2 -hybridized carbon from the graphene lattice to sp^3 -hybridized carbon, giving insight in the degree of functionalization of the nanomaterial. However, it cannot provide any information about the molecule bound to the graphene, making it hard to conclude if the target molecule is attached or not.

Methods

The nanohybrid was characterized by laser desorption/ionization time-of-flight mass spectrometry (LDI-ToF MS) using a Bruker Reflex IV mass spectrometer equipped with a N_2 -laser ($\lambda = 337$ nm). Further, the hybrid was washed with trifluoroacetic acid in DCM:MeOH (1:2, V:V) five times and each of the green filtrates was analyzed via electrospray ionization (ESI) quadrupole ion trap (QIT) MS (esquire 6000, Bruker).

Results

In this study, graphene was covalently functionalized with 5,10,15,20-tetraphenylporphyrin (TPP) via a diazonium salt intermediate. [2] For the first time our investigation employed LDI-ToF MS. LDI revealed the presence of several interesting species on the graphene sample. The molecular-ion peak of the TPP was observed together with other TPP-related signals occurring from side reactions, such as a TPP dimer. Washing the hybrid system several times yielded green filtrates with the color fading from the first to the last washing cycle. ESI MS of the filtrates resulted in signals of TPP together with the byproducts already detected in LDI, decreasing in intensity after each washing cycle. The precipitates from the washing cycles were then investigated again by LDI to check whether TPP is still present on the graphene or not. TPP was still observable in the LDI mass spectrum whereas the signals for the byproducts were significantly reduced, supplying strong evidence that the TPP is indeed covalently attached to the graphene lattice.

Conclusions

LDI-ToF MS was successfully introduced as a new analytic method for the study of covalently functionalized graphene. LDI-ToF MS allowed for the first time the identification of large molecules on a graphene lattice, enabling the detection of products and byproducts, which is not possible with currently applied methods, such as Raman spectroscopy. Furthermore, a combination of ESI-QIT and LDI-ToF MS gave strong evidence that the TPP is indeed covalently attached to the graphene lattice.

Novel Aspect

LDI-ToF MS was introduced as a new analytic method for covalent graphene samples and provided insight into the composition and chemistry of the molecules on the surface.

References

1. Abellán G., Schirowski M., Edenthalhammer K. F., Fickert M., Werbach K., Peterlik H., Hauke F., Hirsch A., J. Am. Chem. Soc., 139, 5175-5182 (2017).
2. Dasler D., Schäfer R. A., Minameyer M. B., Hitzenberger J. F., Hauke F., Drewello T., Hirsch A., J. Am. Chem. Soc., 139, 11760-11765 (2017).

ThP-187 / THE MEASUREMENT OF THE OXIDATION REACTION OF FERROCENES USING LC-MS WITH ESI

Yutaka Okada - Masako Tsuchida

Ritsumeikan University, Department of Applied Chemistry, Kusatsu

Keywords: ESI, LC-MS, ferrocene, oxidation, picrate

Introduction

The oxidation reaction of ferrocenes is one of the typical reactions of ferrocenes. Ferrocene is readily oxidized to the ferrocenium ion under acidic conditions. The oxidation reaction of ferrocene derivatives by air under acidic conditions has been investigated by several researchers. However, the formed ferrocenium ions is not recognized by NMR due to its paramagnetic property. In this paper, LC-MS with ESI was used to detect and quantify the ferrocenium ions.

Methods

The LC-MS was measured using a Shimadzu LCMS-QP8000A. The electric potential was 4.5 kV for positive mode, and -3.0V for negative mode. An OOS column was used under MeOH-H₂O eluent conditions. The oxidation reaction was proceeded by air in 1,2-dichloroethane in the presence of trichloroacetic acid [1]. After the reaction, the solution was washed by water to remove the acid. Picric acid was added to the solution as anion source.

Results

The ESI(+)-MS of ferrocene showed molecular ion peak ($m/z=186$) as base peak, and a fragment ion peak of $m/z=121$. After the oxidation reaction, the mass chromatogram showed two peaks; one is for the unreacted ferrocene, and another is a new peak. The ESI(+)-MS for the new peak showed only molecular ion peak ($m/z=186$). The ESI(-)-MS showed only one peak for picrate anion ($m/z=228$). Therefore, this chromatogram's new peak would be assigned as ferrocenium picrate.

The ESI(+)-MS of ethylferrocene showed molecular ion peak ($m/z=214$) as base peak, and a stable fragment ion peak of $m/z=199$. After the oxidation reaction, the mass chromatogram showed only one peak. The peak was assigned as the oxidized form by ESI(+)-MS and ESI(-)-MS. This is attributed to the high oxidation reactivity of ethylferrocene.

For acetylferrocene, which has low reactivity, small peak was found in the mass chromatogram. The peak was assigned as the oxidized form by ESI(-)-MS.

Conclusions

Based on the LC-MS with ESI, the oxidation reaction of ferrocenes was measured. The oxidized forms were assigned by the ESI(+)-MS, and the picrate ion as the counter anion by ESI(-)-MS. This ESI(-)-MS measurement is useful for the low oxidation reactive ferrocenes such as acetylferrocene.

Novel Aspect

The LC-MS with ESI using picrate ion as the counter anion was an effective method for the measurement of the oxidation reaction of ferrocenes.

References

1. Yamamoto S., Okada Y., Organic communications, 10, 11 (2017).

ThP-188 / MASS SPECTROMETRY-BASED METHODS AS A TOOL FOR MONITORING THE MICRO SCALE CHEMOENZYMATIC PEPTIDE SYNTHESIS ON SOLID SUPPORT

Halina Płóciennik - Piotr Stefanowicz

Keywords: chemo-enzymatic ligation, omniligase, solid phase, LC-MS

Introduction

Synthetic polypeptides and proteins are a particularly useful tools in biochemical research. The synthesis on solid support is limited to only relatively short sequences, while the method of producing recombinant proteins' main restriction lies in incorporating post-translational modifications and unnatural amino acids [1]. A native ligation approach requires cysteine residues in peptide chains and therefore limits the choice of synthetic strategy [2].

Methods

Substrates for the ligation reactions were synthesized on the solid support. Products of synthesis were analyzed using ESI-MS and MS/MS methods. Ligation reactions were performed on the solid support and in solution. Products of ligation were investigated by LC-ESI-MS analysis.

Results

Omniligase is characterized by its unspecific activity. It allows coupling almost without any limitation (except Pro in active center) peptides [3]. Our results show that it is possible to performed ligation on the solid support. One of the most important factor to achieved a product of ligation is chosen appropriate resin. We tested selected solid supports, swellable in water, such as ChemMatrix resin, PEGA resin and TentaGel. The best results were obtained for PEGA resin which was used to the further experiments. Preliminary results indicate the relatively high efficiency of solid phase enzymatic ligation.

Conclusions

In my poster I will compare the results of Omniligase ligation on the solid support and in solution. I will also discuss the factors influencing the efficiency ligation under various conditions.

Novel Aspect

Evaluation of mass spectrometric methods (LC-MS, MS/MS) as a tool for monitoring the micro scale chemoenzymatic peptide synthesis on solid support.

References

1. Guzman F., Barberis S., Illanes A., Electron. J. Biotechnol., 10, 279-314 (2007).
2. Dawson P., Kent S., Annu. Rev. Biochem., 69, 923-960 (2000).
3. Toplak A., Nuijens T., et al., Adv. Synth. Catal., 358, 2140-2147 (2016).

ThP-189 / LC-MS CHARECTERIZATION OF DYNAMIC COMBINATORIAL LIBRARIES OF TASP MOLECULES WITH NOVEL SCAFFOLD

Grzegorz Wolczanski⁽¹⁾ - *Maciej Modzel*⁽²⁾ - *Piotr Stefanowicz*⁽¹⁾ - *Marek Lisowski*⁽¹⁾ - *Thomas J. D. Jørgensen*⁽²⁾

Faculty of Chemistry, University of Wrocław, Joliot-Curie 14, 50-383 Wrocław, Poland (1) - University of Southern Denmark, Department of Biochemistry and Molecular Biology, Campusvej 55, 5230 Odense Denmark (2)

Keywords: TASP, Dynamic combinatorial libraries, LC-MS, synthetic helical proteins

Introduction

TASP molecules are artificial protein models formed by attaching peptide chains to a cyclic scaffold [1]. Recently we reported the preparation of DCLs of helical proteins models basing of aromatic disulfide

chemistry [2]. We present preparation of TASP dynamic combinatorial libraries (DCLs) basing on a new scaffolds, which are created by formation of disulfide bonds between *N*-(2,6-dimercapto-1,3,5-triazin-4-yl)peptides. DCLs were analyzed by LC-MS.

Methods

We synthesized our peptide conjugates using the standard SPPS protocol. Libraries of TASP molecules have been created from *N*-(2,6-dimercapto-1,3,5-triazin-4-yl)peptides (or their dimers with a disulfide bridge) *via* oxidation by atmospheric oxygen or CuCl₂ in 2-propanol. Selective formation of a cystine residue has been performed by adding I₂. We used LC-MS and IM-MS to characterize TASP libraries. Identity of the isolated compounds has been confirmed by HR-MS.

Results

Oxidation of *N*-(2,6-dimercapto-1,3,5-triazin-4-yl)peptides by atmospheric oxygen in aqueous solution does not occur. The peptides incubated in a more hydrophobic environment of 2-propanol produced TASP molecules by formation of disulfide bridges between 2,6-dimercapto-1,3,5-triazine moieties. Subsequent oxidation of the peptide conjugates containing an *N*-terminal Cys residue with an attached dimercapto triazine moiety allows controlling the size of the macrocyclic scaffold. Similarly to systems proposed by Otto *et al.* [3] and our previous results [2], the peptide conjugates gave mainly trimers and tetramers. Formation of disulfide bond between Cys residues attached to triazine ring affects the size of oligomers. Used peptide sequence were designed as short amphipathic α -helices. Dynamic combinatorial libraries of the oligomers were analyzed by nLC-MS and TWIMS-MS. Analysing of chromatographic behavior and experimental arrival times, we concluded that tetramers may occur as two main topoisomers - C₄ and C_i symmetry.

Conclusions

We developed a new type of scaffold based on triazine moieties connected by disulfide bonds into cyclooligomers. Our method is useful in designing dynamic combinatorial libraries of TASP. The libraries of synthetic proteins may be simply characterized using LC-MS and IM-MS. Fixing of the relative positions of the two triazine rings by subsequent formation of disulfide bonds between attached Cys residues affects the size of the macrocyclic scaffold.

Novel Aspect

We developed first dynamic combinatorial libraries of TASP molecules with scaffolds based on oligo triazine disulfides, which were characterized by LC-MS and TWIMS-MS.

References

1. Tuchscherer G., Mutter M., 41, 197 (1995)
2. Wołczański G., Cal M., Waliczek M., Lisowski M., Stefanowicz P., Chemistry, 2018, <https://doi.org/10.1002/chem.201800187>
3. Komáromy D., Stuart M. C. A., Santiago G. M., Tezcan M., Krasnikov V. V, Otto S., J. Am. Chem. Soc., 139, 6234 (2017)

ThP-190 / ELEMENTAL ANALYSIS BY AU-NANOPARTICLES ENHANCED LA-ICP-MS

Fabrizio Mastroiocco⁽¹⁾ - **Lorena Carla Giannossa**⁽²⁾ - **Annarosa Mangone**⁽¹⁾

Università di Bari, Department of Chemistry, Bari, Italy (1) - Università di Bari - Department of Economics, Bari, Italy (2)

Keywords: Laser Ablation, Nanoparticles, LA-ICP-MS

Introduction

A preliminary study of nanoparticles-enhanced LA-ICP-MS is presented. It represents an analytical strategy

to improve the analytical performance of the LA-ICP-MS in terms of sensitivity and LOD, preserving its excellent characteristics (simplicity, rapidity, etc.), without any changes in the experimental set-up.

Methods

Some drops of Gold Nanoparticles (AuNPs) colloidal dispersion were deposited on the sample's surface, the solvent evaporated before the analysis. Spherical AuNPs dispersions in aqueous citrate buffer (Sigma Aldrich Co.) or naked PLAL (Pulse Laser Ablation in Liquid) [1] AuNPs of different size were used on certified standard target (Copper, Titanium, Bronze, Brass, Silicon and Glass).

Results

A considerable increase of the measured signal intensity was observed in the presence of AuNPs. By analogy with results already obtained with Laser Induced Breakdown Spectroscopy [2,3], when a critical number of NPs are deposited on the target surface, the laser pulse electromagnetic field induces the collective oscillation of the conduction electrons of the NPs that in turn results in a strong enhancement of the field. The improvement is probably induced by the best quality of the nanoparticles' agglomerate generated during the laser ablation.

Metallic elements show enhancement also in non-conductive matrices. Different metallic elements show different enhancement in the same matrix, as well as the same element shows different enhancement in different matrices. It should be noted the AuNPs can be completely removed after the analysis.

Conclusions

The existence of an improvement of sensitivity in LA-ICP-MS signal via AuNPs allows to identify and quantify elements without damaging significantly the sample's surface unlike conventional LA-ICP-MS.

Novel Aspect

The AuNPs allow to improve the sensitivity of LA-ICP-MS technique almost in a nondestructive mode, which can be extremely useful in many fields as forensic sciences and cultural heritage.

References

1. De Giacomo A., Dell'Aglio M., Santagata A., Gaudio R., De Pascale O., Wagener P., Messina G. C., Compagnini G., Barcikowski S., Phys. Chem. Chem. Phys. 15, 3083–3092 (2013).
2. De Giacomo A., Gaudio R., Koral C., Dell'Aglio M., De Pascale O., Analytical Chemistry, 85, 10180-10187 (2013).
3. De Giacomo A., Gaudio R., Koral C., Dell'Aglio M., De Pascale O., Spectrochimica Acta Part B 98, 19–27 (2014).

ThPS-S10 - Novel food; nanoparticles in food and beverage; future applications

Chair: Stefan Weigel

ThP-191 / A NEW BY-PRODUCT ENRICHED IN PHENOLIC COMPOUNDS, EVALUATED BY UHPLC-HRMS, FROM OLIVE MILL WASTEWATER

Giuseppe Avellone⁽¹⁾ - ***Vita Di Stefano***⁽¹⁾ - ***Maria Eugenia Novara***⁽²⁾ - ***Rosa Pitonzo***⁽²⁾

UDS Palermo, STEBICEF, Palermo (1) - UDS Palermo, ATeN Center, Palermo (2)

Keywords: Olive mill wastewater, by- product, UHPLC-HRMS, hydroxytyrosol

Introduction

Olive mill wastewater (OMW), a by-product obtained during olive oil extraction, represents an important environmental problem for the high toxicity, the resistance to biodegradation and the high volumes processed. Instead of the pollution, OMW is a rich source of bioactive compounds, polyphenols, which could be found up to 98% of the total content of olive fruits [1]. Hydroxytyrosol is the main detected phenolic alcohol [2, 3], widely known due to the natural antioxidant and anti-inflammatory properties [4]. The aim of the present study is to explore and compare the qualitative and quantitative phenolic profile of two OMW treated by two different techniques: freeze-drying and spray-drying.

Methods

Two samples of OMW, Ava and Bava, from two different Sicilian oil mills and generated by three phase continuous systems, were centrifuged, the supernatant was dismissed while the solution was subjected to freeze-drying (FD) and spray drying (SD). The experimental conditions of these techniques were optimized. The phenolic compounds were extracted by ethyl acetate and, finally, were analyzed by ultra-high-performance liquid chromatography-heated electrospray-high resolution mass spectrometry (UHPLC-HESI-HRMS).

Results

Qualitative and quantitative phenols composition were evaluated by UHPLC-HESI-HRMS. For freeze-drying and spray drying samples were identified and quantified about eleven polyphenols. This new HRMS analyzer allows to characterize in complex mixtures small molecules with confident and rapid monitoring. The qualitative and quantitative analysis showed an important complexity of this wastes, in which, thanks to the new HRMS analyzer, it was able to determinate the high concentration of polyphenols belonging to phenyl acids, phenyl alcohols and secoiridoids derivatives. Hydroxytyrosol is the most important phenolic compounds, extracted in FD and SD samples, and the concentration, up to 18,11 mg/g (in FD OMW). This amount was determinate for the first time with this high concentration. Elenolic acid and hydroxy-O-decarboxymethyl-oleuropeina aglycone were the second phenolic compounds most represented: 6,57 mg/g in Bava FD and 1,87 mg/g in Bava SD and the second one was 4,48 mg/g in Bava FD and 1,82 mg/g in Bava SD. Freeze-dry was the best methods to enrich in phenolic composition the OMW samples especially in hydroxytyrosol.

Conclusions

The possibility of using and recycling bioactive compounds from agroindustrial wastes would reduce pollution caused by the latter, and would also allow to realize pharmaceutical, cosmetic, food and feed products enriched with functional molecules. In this view, new technologies, proposed to treating of OMW, have been used to valorise these wastes as a low-cost starting material as predict by new bio economy. The recovered polyphenols obtained by freeze-drying could lead towards better usage of by-products from olive oil production.

Novel Aspect

These techniques allow to recycle the OMW, only, through physic treatment without any chemicals extraction or processes. For SD samples was added maltodextrin, widely known as food additive, which necessary to maintain the stability over the time.

Therefore, the agricultural industries could implement the FD or SD systems in order to obtain valuable health-products while containing costs, especially for SD.

References

1. Rodis P. S., Karathanos V. T., Mantzavinou A., Partitioning of olive oil antioxidants between oil and water phases. *J. Agric. Food Chem*, 50, 596-601 (2002).
2. DellaGreca M, Previtiera L, Temussi F, Zarrellia, Low molecular-weight components of olive oil mill wastewaters. *Phytochem Anal* 15:184–188 (2004).
3. H. K. Obied, P. D. Prenzler, K. Robards, Potent antioxidant biophenols from olive mill waste. *Food*

Chemistry 111: 171-178, (2008)

4. A. Cardinali, N. Cicco, V. Linsalata, F. Minervini, S. Pati, M. Peralice, N. Tursi, V. Lattanzio, Biological Activity of High Molecular Weight Phenolics from Olive Mill Wastewater. *J. Agric. Food Chem.*, 58, 8585–8590 (2010).
5. M. D. Di Mauro, R. C. Giardina, G. Fava, E. F. Mirabella, R. Acquaviva, M. Renis, N. D'Antona. Polyphenolic profile and antioxidant activity of olive mill wastewater from two Sicilian olive cultivars: Cerasuola and Nocellaraetnea. *Eur Food Res Technol*, 243, 1895–1903 (2017)

ThP-192 / OPTIMIZATION OF ESPRESSO COFFEE EXTRACTION WITH DIFFERENT PARTICLE SIZE DISTRIBUTION AND ANALYSIS THROUGH GC-MS AND HPLC-VWD

Gulzhan Khamitova^(1,2) - ***Simone Angeloni***^(1,2) - ***Giovanni Caprioli***⁽¹⁾ - ***Gianni Sagratini***⁽¹⁾ - ***Sauro Vittori***^(1,2)

University of Camerino, School of Pharmacy, Camerino (1) - International Hub for Coffee Research and Innovation, University of Camerino, Via S. Agostino 1, 62032 Camerino, Italy (2)

Keywords: Espresso Coffee (EC), Particle Size Distribution (PSD)

Introduction

EC extraction depends deeply on ground and roasted coffee particles produced by the grinding and brewing process. To extract tempting aroma and taste in EC, the grinding process is a crucial step. The taste and flavor change in important manner owing to the particle size of ground coffee [1]. This research aims at optimizing the impact of different particle size on aroma and bioactive compounds.

Methods

Quantitative and qualitative analyses were carried out on Arabica and Robusta cultivars. Different particle size average (200 and 400 microns) were used for extraction. Volatiles and bioactive compounds analyses were carried out through GC-MS and HPLC-VWD, respectively. The usual espresso machine settings (9 bars and 92°C) were kept [2].

Results

EC extraction times, kinetics of extraction and its dependence on water pressure and temperature were controlled during the analysis. Volatile compounds, determined in different particles, were divided into family groups and their characteristics were identified [3]. According to literature[4], particle size influences the amount of extracted bioactive compounds in EC. To study the effect of the particle size, EC samples were processed analyzing caffeine, nicotinic acid and trigonelline. EC extraction were in triplicate for each condition for each particle size. results confirmed a good extraction efficiency of caffeine, which accounted for 170 mg and 90 mg per cup, on 200 and 400 microns, respectively. Hence, amount of bioactive compounds were increased extracting smaller particle sizes.

Conclusion

EC extraction of different particle sizes were analyzed for their content in bioactive components and volatiles normally found in EC. Extraction optimization will be developed by modifying of extraction variables and further studies will be conducted.

Novel Aspect

The feasibility study of particle size effect on extraction of EC has not been well established yet. Understanding particle size impact on EC will help in setting optimal extraction conditions.

References

1. B. Folmer, I. Blank, A. Farah, P. Giuliano, D. Sandres, C. Wille, *The Craft and Science of Coffee*, London,

UK: Elsevier (2017).

2. G. Caprioli, M. Cortese, G. Cristalli, F. Maggi, L. Odello, M. Ricciutelli, G. Sagratini, V. Sirocchi, G. Tomassoni, S. Vittori, Optimization of espresso machine parameters through the analysis of coffee odorants by HS-SPME-GC/MS, *Food Chemistry*, 135, 1127-1133 (2012).
3. J. A. S. Lopez, M. Wellinger, A. N. Gloess, R. Zimmermann, C. Yeretian, Extraction kinetics of coffee aroma compounds using a semi-automatic machine: On-line analysis by PTR-ToF-MS, *International Journal of Mass Spectrometry*, 401, 22-30 (2016).
4. M. Kuhn, S. Lang, F. Bezold, M. Minceva, H. Briesen, Time-Resolved Extraction of Caffeine and Trigonelline from Finely-Ground Espresso Coffee with Varying Particle Sizes and Tamping Pressures," *Journal of Food Engineering*, 206, 37-47 (2017).

ThP-193 / STUDYING OF PHENOTYPIC AND PHYSIOLOGICAL INFLUENCES OF ABIOTIC STRESSORS (TiO₂ AND EXCESS IRRADIATION) ON TOMATO PLANTS

Jung Ah Ko

Korea Institute of Toxicology, Future Environment Research Center, Jinju, Republic of Korea

Keywords: environmental pollutant, Titanium dioxide (TiO₂), phytotoxicity, climate change, nanoparticle behavior

Introduction

Titanium dioxide (TiO₂) nanoparticles are widely used for skin care, food additives, solar cells, and photocatalysis. However, high demand of TiO₂ raise chance of the releasing into the environment as a pollutant. [1] Since nanoparticles including TiO₂ can arise toxicological effect based on their physicochemical activity, [2-3] studying and expecting the possible influence to the human health is important issue in recently.

Methods

Herein, we studied about influences of TiO₂ when it exposed to tomato plants as a pollutant under excess light condition, for simulating behaviors of TiO₂ under climate change. Plant cultivation was performed in hydroponic system, and phenotypic and physiological analysis were performed to monitor and to suspect any possible influences on the human health.

Results

Accumulation of TiO₂ in leaves were analyzed by qualitative and quantitative measurements with STEM, XRD, and ICP-MS. Through the obtained result, we have confirmed clear evidence of TiO₂ presence based on the plant parts, *i.e.*, high accumulation of TiO₂ was detected in leaves and root, while no storage was analyzed in stem.

Other physiological effects on the model plant were exhibited, for example, increased level of oxidative enzyme activity and antioxidant level, color and thickness change of leaves, enhancement on photosynthetic rate, and reproductivity, with respect to the abiotic stressors of TiO₂ and excess irradiation. Presence of TiO₂, especially, lead to increase of the photosynthetic efficiency by improving the activity that related with photosynthesis, such as quantum yield of PSII, electron transport rate, stomatal conductance and so far.

Conclusions

Abiotic stressors of TiO₂ and excess light induced more oxidative stress in comparison with control, and the increased stress expected arisen by direct physical damage as well as enhancement of photosynthetic efficiency. However, we confirmed raised antioxidant level, and higher growth efficiency and fruits production, of which can improve quality of fruit and fruit yield in agricultural field from the tomato plants grown under the stressful environment.

Novel Aspect

This study focused on not only phenotypic influence but also physiological effect of tomato plants to expect possible effect on human health.

References

1. Arturo A. K., Anastasiya L., Environ. Sci. Technol. Lett., 1, 65 (2014)
2. Bradley J. C., Raven B., Courtney K., J. Nanopart. Res., 14, 913 (2015)
3. Camille L., Julien L., Nathalie H. -B., Hicham K., Barbara F., Anne-Marie F., François B., Marie C., Sci. Total Environ., 431, 197 (2012)

ThPS-S11 - Informatics tools and data analysis

Chairs: Peter O'Connor, Giuseppe Avellone

ThP-194 / TAMI SOFTWARE FOR THE PROVISION OF ELEMENTAL FORMULA FROM YOUR SINGLE QUADRUPOLE GC-MS OR LC-MS DATA

Aviv Amirav - Alexander B. Fialkov - Tal Alon

School of Chemistry, Tel Aviv University, Tel Aviv 69978, Israel

Keywords: Elemental formula from quadrupole MS

Introduction

In order to significantly improve sample identification by GC-MS we developed the Tal-Aviv Molecule Identifier (TAMI) software. TAMI inverts the molecular ion isotopomeric pattern into elemental formulae and automatically confirms or rejects NIST library identification for improved confidence in the identification. When NIST identification fails, TAMI provides the sample compound elemental formula from standard centroid data files.

Methods

When the quadrupole MS scans in profile mode, TAMI improves the mass accuracy (typically <50 ppm) and combined with isotope distribution analysis it usually provides the correct elemental formula as number one in its hit list.

The TAMI software was extensively explored with Agilent 5975/5977 GC-MS, Aviv Analytical 5977-SMB GC-MS with Cold EI that provides enhanced molecular ions, Thermo ISQ Chromeleon and Advion CMS LC-MS data files with the following main benefits:

Preliminary data

A) Automatic linking with the library. TAMI automatically confirms or rejects the NIST library identification via independent isotope abundance analysis.

B) Independent provision of elemental formula from standard centroid data files. In case of rejection of the library identification such as when the compound is not in the library, TAMI independently provides a list of elemental formulae with declining order of matching to the experimental data.

C) Simple to use. The linking of TAMI with the NIST library is automated while the provision of elemental formula requires only one mouse click.

D) Improved mass accuracy. When the quadrupole MS scans in profile mode TAMI improves the quadrupole mass accuracy via the use of special weighted average algorithm and mass calibration with PFTBA to be typically better than 50 ppm.

E) Special compatibility with Thermo Chromeleon Data Files. The Thermo ISQ in combination with Chromeleon software provides <50 ppm mass accuracy using the TAMI PFTBA based mass tune function with standard centroid data files.

Novel Aspect

Software for the provision of elemental formula from unit resolution quadrupole GC-MS or LC-MS based on isotope abundance analysis and improved mass accuracy

ThP-195 / THE DATABASE FOR CROSS-PLATFORM MINING BIOMARKERS IN THE BRAIN TUMOR TISSUE

Konstantin Bocharov⁽¹⁾ - *Anatoly Sorokin*⁽¹⁾ - *Vsevolod Shurkhay*⁽²⁾ - *Igor Popov*⁽¹⁾ - *Evgeny Zhvansky*⁽¹⁾ - *Alexander Potapov*⁽²⁾ - *Evgeny Nikolaev*⁽³⁾

Moscow Institute of Physics and Technology, Laboratory of Ionic and Molecular Physics, Dolgoprudny (1) - Burdenko Scientific Research Neurosurgery Institute, Moscow (2) - Skolkovo institute of science and technology, Skolkovo (3)

Keywords: mass-spectrometry, database, biomarkers, brain tumors, high and low resolution

Introduction

One of the most challenging problems when using highly specialized mass spectrometry instruments and ionization techniques is both the unification of the raw data obtained and the interpretation of the final result of their processing. So, we present here an approach using database for comparing the high and low resolution spectra obtained by various ionization methods on different mass spectrometers and intended for identification of different tumor types.

Methods

Tumor tissues were dissected during neurosurgery. The spectra were measured by: a) needle electrospray ionization [1] on ICR-based HR mass spectrometer LTQ FT Ultra; b) direct spray of lipids extracted from tissue by Folch method on Thermo LTQ-Orbitrap XL and Bruker Esquire 3000 plus. The high and the resolution was 56,000 and 1,000 at 800 Th, respectively. Information about ionization method, instrument type, tumor type, etc. was stored in NoSQL database.

Results

A database of over 300 samples of different tumor tissues and surrounding tissues from over 250 patients was collected in the study. All of the samples were provided by the Burdenko Neurosurgical Institute. Collected spectra were filtered and stored in column-oriented database MonetDB. User interface allows for visualisation, comparison, transformation (PCA, ICA etc.) of the spectra. The integration with R/Bioconductor makes whole spectra of statistical analysis and machine learning techniques readily available from API as well as from user interface. To demonstrate the infrastructure performance 9 classifiers were trained on 20 astrocytoma and 9 unmodified brain samples with different settings. Comparative analysis of the results obtained from spray-from-tissue ionization experiments and extracts analysis shows that both methods demonstrate similar series of peaks. All major peaks are found in both types of spectra. Peaks in low resolution usually are larger when normalized to TIC.

Conclusions

Vast majority of selected features distinguishing the tumor tissue from the adjacent tissue was found in both high- and low-resolution spectra. We have created mapping schema, which allows using classifier trained to separate tumor samples from non-tumor brain samples on high-resolution spectra with data from low-resolution spectra with a minor decrease in accuracy.

Novel Aspect

Our cross-platform approach is a huge step to the unification of the data obtained at various types of mass spectrometry platforms, in particular, for LC-MS.

References

1. Kononikhin A. et al., Anal Bioanal Chem, 407, 7797-805 (2015).

Acknowledgement

Supported by the RSF (16-15-10431).

ThP-196 / COMBINING GC-MS WITH MACHINE LEARNING FOR NON-INVASIVE CANCER AND DISEASE DIAGNOSTICS

Michael Cauchi

University of Limerick, Mathematics & Statistics, Limerick, Ireland

Keywords: Metabolomics; GC-MS; Machine learning; Biomarker discovery; Medical diagnostics

Introduction

The use of machine learning in biomedical applications has grown. Gold standard techniques are invasive, costly and time-consuming. A non-invasive approach measures the volatile analytes present in patients' urine or faeces via gas chromatography mass spectrometry (GC-MS) in conjunction with optimised machine learning models that can distinguish between ill and healthy patients. Bladder cancer [1] and Crohn's disease [2] are presented.

Methods

Patients diagnosed with bladder cancer or Crohn's disease were recruited and donated urine or faecal samples respectively. Samples were processed and passed through a GC-MS instrument. Each sample generated a data matrix from which the total ion count (TIC) chromatogram was attained. Normalisation was followed by peak alignment via correlation optimized warping (COW). The data were subjected to machine learning algorithms for diagnosis and biomarker discovery.

Results

The three machine learning algorithms employed were partial least squares discriminant analysis (PLS-DA), support vector machines (SVMs) and random forests (RFs). The models were optimized via a bootstrapping process in order to ensure confidence in the results. The bladder cancer data (from urine) consisted of three control groups (no disease, infection and urological disease) and three cancer grade groups. Of great interest was to see if urological diseases could be distinguished from the cancer grades: 77%, 79% and 74% overall accuracies were attained for PLS-DA, SVMs and RFs respectively. The other control group comparisons attained higher accuracies (*i.e.* >85%). The Crohn's disease data (from faeces) also produced a very good overall accuracy of 85%. Interrogation of the optimal PLS-DA models for the bladder cancer and Crohn's disease suggested significant retention time peaks which led to the identification of possible markers.

Conclusions

The application of GC-MS in conjunction with machine learning techniques offers a highly feasible, non-invasive procedure for cancer and disease diagnosis. In the case of bladder cancer there was an improvement on the gold standard technique of urine cytology. In the case of Crohn's disease, results compared favourably with the gold standard technique of colonoscopy.

Novel Aspect

Combining GC-MS with machine learning offers a non-invasive approach to cancer and disease diagnostics. Applications to other diseases/instruments are highly feasible.

References

1. Cauchi M, et al., Analytical Methods, 8, 4037-4046 (2016).
2. Cauchi M, et al., Metabolomics, 10, 1113-1120 (2014).

ThP-198 / INTEROPERABLE AND SCALABLE METABOLOMICS DATA ANALYSIS WITH MICROSERVICES

Payam Emami Khoonsari⁽¹⁾ - ***Stephanie Herman***⁽¹⁾ - ***Marco Capuccini***^(2,3) - ***Matteo Carone***⁽³⁾ - ***Anders Larsson***⁽⁴⁾ - ***Joachim Burman***⁽⁵⁾ - ***Steffen Neumann***⁽⁶⁾ - ***Christoph Steinbeck***⁽⁷⁾ - ***Kim Kultima***⁽¹⁾ - ***Ola Spjuth***⁽³⁾ ***and the PhenoMeNal project partners***

Department of Medical Sciences, Clinical Chemistry, Uppsala University, Uppsala (1) - Department of Information Technology, Uppsala University, Uppsala (2) - Department of Pharmaceutical Biosciences, Uppsala University, Uppsala (3) - National Bioinformatics Infrastructure Sweden, Uppsala University, Uppsala (4) - Department of Neuroscience, Uppsala University, Uppsala (5) - Department of Stress & Developmental Biology, Leibniz Inst. of Plant Biochemistry, Halle (6) - Friedrich-Schiller-University, Jena (7)

Life scientists working with big data such as genomics, proteomics and metabolomics are beginning to experience numerous difficulties of handling, processing and moving information that has not been experienced before. These intensive computational tasks cannot be tackled by regular laptops and workstations. The high demand for more resources and increased flexibility has even turned the well-known shared high-performance computing unresponsive.

Cloud computing together with microservices architecture offer a compelling alternative to these old models, with the possibility to instantiate and configure on-demand resources such as virtual computers, networks, and storage, together with operating systems and isolated software tools.

We have developed an architecture which uses components for data analysis encapsulated as microservices connected into operating computational workflows. This solution provides users with complete, ready-to-run, reproducible and scalable data analysis environment that can be easily deployed on desktop computers as well as public and private clouds, without requiring special IT skills from the user.

We showcase the capability of this architecture using four demonstrators: in the first demonstrator we reproduce a large metabolomics data analysis and show how it can be effortlessly scaled up on cloud computing resources with little efficiency loss. In this case, 1092 computationally intensive tasks were performed in less than four hours which otherwise could have taken days to finish on a workstation. In the second demonstrator we provide new mass-spectrometry data from multiple sclerosis patients, and present for the first time, to the best of our knowledge, a complete start-to-end analysis of untargeted mass-spectrometry and identified novel biomarkers. The last two demonstrators describe NMR workflow and fluxomics workflows, emphasizing that microservice architecture is domain-agnostic.

The four demonstrators show the versatility, applicability and scalability of our method in metabolomics; however, the methodology can be applied to all scientific disciplines, paving the way towards for large-scale integrative science.

ThP-199 / VALIDATION OF ANALYTICAL METHODS FOR THE DETERMINATION OF RESIDUES IN FOOD: THE STUMBLING BLOCK OF DATA PROCESSING

Roberta Galarini

Istituto Zooprofilattico Sperimentale dell'Umbria e delle Marche "Togo Rosati", Center for Method Development and Validation, Perugia

Keywords: method validation, data processing, residues, food control

Introduction

Starting from the late 90s, the ISO accreditation standard along with the EU rules required a great effort to validate methods and to process the relevant data in order to demonstrate the fitness for purpose of internal developed procedures. Today, this aspect is even more critical because, thanks to MS detectors, multiclass methods including dozens of analytes are currently applied.

Methods

The free software package called ADVerSE (Analytical methods Validation results Speedy Elaboration) has been developed in R environment [1]. ADVerSE produces the Validation Report processing the validation data put into a predefined Excel file. In the same file the legal requirements fixed for each coupling analyte/matrix (maximum permitted levels) are inserted, too.

Results

The first part of Validation Report (.doc) generated by the ADVerSE package includes the descriptive statistics of results, *i.e.* mean recovery and precision at each validation level for each analyte. Linearity (linear regression) is assessed applying Mendel and lack-of-fit tests. In the second part, the precision is studied as a function of concentration; limits of detection and quantification are estimated using the results obtained at the lowest concentration of the validation study. With regards to the assessment of compliance with legal limits, the measurement uncertainty is estimated applying the top-down approach [2]. Finally, according to the specific requirements stated by Commission Decision 2002/657/EC (veterinary drug residues in food), the limit of decision (CC α) and detection capability (CC β) are calculated both for banned and permitted substances [3].

Conclusions

The created Report allowed the overcoming of the long data processing saving time and minimizing calculation mistakes. The few seconds required for its generation favor an immediate evaluation of method performance characteristics helping the analyst to eventually improve her/his analytical procedure or the experimental design adopted for the validation study.

Novel Aspect

This is the first package able to produce a Validation Report usable in the various legislative contexts of the official control of residues in food.

References

1. The Comprehensive R Archive Network. <https://cran.r-project.org>
2. EURACHEM/CITAC. Guide Quantifying Uncertainty in Analytical Measurement. Third Edition (2012)
3. Commission Decision 2002/657/EC. Official Journal of the European Communities. L221, 8-36 (2002)

ThP-200 / INTERFACING THIRD PARTY SOFTWARE APPLICATIONS TO MASS SPECTROMETRY DATA SYSTEMS

Jeff Goshawk - Mike McCullagh - Russell Mortishire-Smith

Waters, Scientific Operations, Wilmslow

Keywords: Application Programming Interface, Library Generation

Introduction

The desire to extract mass spectrometry data from commercial systems is becoming increasingly popular.

Customers of such systems wish to use data in bespoke applications and incorporate data into existing repositories. An Application Programming Interface (API) is proposed that allows access to raw and processed mass spectrometry data contained within the instrument platform. An external application using the API to create a library of compounds is discussed.

Methods

A representational state transfer (REST) API uses HTTP requests to move, change, create and delete data. Data acquired on an ion-mobility enabled instrument for a selection of 500 natural products were processed in the UNIFI mass spectrometry data system. An application, written using the UNIFI API, was employed to automatically interrogate the processed data, then to aggregate the reduced outcomes, to produce a library of compounds.

Results

Manual analysis and collation of mass spectrometry data to create a library of compound information is a laborious task that is prone to error. An application, external to the UNIFI mass spectrometry data system, has been written that accesses data through the UNIFI API and automatically extracts processed data to produce compound library content. The library content produced by the application is rich in content and can optionally consist of retention times, multiple adducts, together with their individual collisional cross section measurements, and fragment ion m/z values which have been generated using an in silico structure fragmentation tool and matched to measured high energy ions. All quantities derived can be averaged over multiple replicate injections. The application transforms a manual process, which often takes weeks or even months, and is open to transcription errors, into an automated process that can be performed error-free in a matter of hours.

Conclusions

The UNIFI Application Programming Interface allows direct access to raw and processed data within UNIFI. A compound library generation application has been presented which accesses data through the UNIFI API and then uses the data to derive accurate library content. A library of commercially natural product compounds has been created and contains information on multiple adducts, collisional cross sections and fragment ions.

Novel Aspect

Direct access to mass spectrometry data for use with third party applications. Automatic production of mass spectral library content through the analysis of processed data.

ThP-201 / METABOLOMIC ANALYSIS OF HIGH DENSITY MS DATA

Andrea Luterova - Tomas Furst - Miroslav Strnad - Jiri Gruz

Palacky University, LGR, Olomouc

Keywords: LC-MS, metabolomic analysis

Introduction

Metabolomic analysis allows to systematically identify and quantify metabolites from a biological sample. Non-targeted analysis which refers to the analysis of all metabolites is typically performed by liquid chromatography-mass spectrometry (LC-MS). The resulting 3D-chromatogram contains thousands of features that need to be properly analysed and processed. This can only be achieved by automated computational approaches.

Methods

MATLAB, a multi-paradigm numerical computing environment (MathWorks, Inc.).

Results

A software for automated data evaluation was developed in our laboratory by using Matlab programming language. Besides basic chromatogram alignment, various post-processes were designed and developed. One of them is grouping of so called “friendly ions” that consist of isotopic, adduct and fragment ions (originating from a single metabolite). Another post-processing approach is the mass correction which is based on known background ions or identified metabolites. This allows lockmass-like correction with no measured lockmass data. The third post-processing method is an algorithm for automated molecular peak identification. This approach is based on identification of typical fragments and adducts under defined experimental conditions. As it relies on correct determination of friendly ions, its effectivity is also affected by retention time stability and separation effectivity.

Conclusions

The developed algorithms can help to process high density metabolomic data with improved accuracy and effectivity.

Novel Aspect

The described algorithms take advantage of the access to raw MS data which greatly improves their effectivity when compared to current methods.

Acknowledgement

This work was supported by the Czech Science Foundation (project GA17-06613S).

ThP-202 / A STRATEGY TO ANALYZE THE SPECTRAL DATASETS OF UNKNOWN COMPOUNDS USING MASS SPECTRAL SIMILARITY NETWORK

Eisuke Hayakawa

Okinawa Institute of Science and Technology, Okinawa Institute of Science and Technology, Okinawa

Keywords: Informatics, Metabolomics, Natural Product, Fragmentation, Structure Analysis

Introduction

In mass spectrometry-based small compound analysis, fragment spectra have been widely used to deduce the structures of analyte compounds. In recent years, approaches such as reference spectrum library searching and structure database-assisted spectrum matching have been used in many applications. However, systematic and high-throughput structure annotation for large scale spectral datasets is still quite difficult to accomplish.

Methods

Mass spectral similarity network, also called as MS/MS spectral similarity network or molecular networking, is a method to assess and visualize the similarities of fragment spectra of structurally related compounds. Since similarities of fragment spectra can be correlated to the structural similarities of the corresponding compounds, this network is a useful tool to detect structural similarity of analyte compounds.

Results

We have developed a data analysis framework employing mass spectral similarity network to organize, classify, (sub)structurally annotate and visualize the complex fragment spectral dataset of unknown compounds. Reference fragment spectra is acquired from public spectral library such as MassBank and used to create “reference” spectral similarity network. Further, in-silico fragment spectra are generated by applying fragmentation models to the subsets of chemical structure database to create “in-silico” mass spectral similarity networks. These spectral similarity networks are further divided into user-defined classes such as chemical classes, substructures, source types. The multi-layer representation of spectral similarity

networks effectively visualizes and categorizes the spectra of unknown compounds in samples. Further, network informatics and cheminformatics methods are applied the network structure to extract and assign the structural information, providing the structural information of the unknown compounds.

Conclusions

This versatile method can be applied to mass spectral datasets containing a wide range of small molecules. We demonstrate the application of this method using case studies.

Novel Aspect

A novel mass spectral similarity-based method to detect, classify and annotate the structures of unknown compounds in large-scale datasets of fragment spectra.

ThP-203 / A LOW PULSE CURRENT MEASUREMENT CIRCUIT FOR CURRENT MONITOR SYSTEM AT MASS SPECTROMETER

ZiXing He - Lihua Zhai - ChaoHui He - QingMin Zhang - Zhiming Li - Xianglong Yuan

Xi'an Jiaotong University, Xi'an Jiaotong University, Xi'an

Keywords: reverse charging; constant current source; blast pulse; mass spectrometer

Introduction

A pulse beams at 105 particles per second have been generated by laser Resonant ionization at 10KHz. Beam current is normally measured with the aid of Faraday cup. However, the signal is difficult to be monitored by Faraday cup amplifier, which was developed to monitor the beam current at pA range. Simultaneously, a discriminator is inappropriate, due to the time of death [1] and response rates. Therefore, What has been designed is a low pulse current measurement circuit, which is used for current monitor system at mass spectrometer.

Methods

A low pulse current measurement circuit, applying synchronous charging module and reversed charging constant current source module, was designed for current monitor system at mass spectrometer. This circuit turns the strength of measuring signal into the amplitude of voltage [2]. As a modified charge integration circuit for periodic blast pulse, a reverse charging circuit is designed [3], aiming to discharge integral capacitor for the avoidance of saturation. At the same time, it requires to improve the efficiency of measurement. So, a controllable constant current source model will be designed.

Results

In circuit simulation once the integral capacitor's voltage exceeds given threshold of 12V, the comparator will generate a pulse signal for discharging. Then the integral capacitor's voltage will drop rapidly. At the end of discharging process, the integral capacitor's voltage just drops to -10V. At the same time, the discharged signal is generated by the comparator as a counting signal sent to the counter. In experiments, A modified charge integration circuit for periodic blast pulse is verified by measuring the synchronous pulse signals that are generated by laser ionization. Verification results show that the experiment error is 2%.

Conclusions

By calculations, simulations and experiments, The results show that a charge measurement circuit for periodic blast pulse is very effective. It is also shown that the count rate of periodic blast pulse signals between 103 and 107 periodic blast pulse can be measured by this circuit. At the same time, the circuit is effective in the continuous measurement of the periodic blast pulse signal. Measurement precision of this circuit is mainly affected by the slew rate of operational amplifier and the bypass capacitors of the input

signal [4]. So, it is feasible to further improve the accuracy of measurement by perfecting the circuit.

Novel Aspect

A modified charge integration circuit for periodic blast pulse, a reverse charging circuit is designed and dynamic range of the ion monitor is effectively expanded.

References

1. Chen Yu, Ma Cheng-xing, High count-rate single channel pulse height analyzer. University of Science and Technology of China, 20, 10-12 (in China) (2000).
2. Alam, A. H. M. Z., et al., Design of Capacitance to Voltage Converter for Capacitive Sensor Transducer. American Journal of Applied Sciences, 7, 1353-1357 (2010).
3. Anania, M. P., et al., Design of a plasma discharge circuit for particle wakefield acceleration. Nuclear Instruments & methods in Physics Research A, 740, 193-196 (2014).
4. Thakker, R. A., et al., Operational amplifier having improved slew rate, US Patent 8,089,314 B2 (2012).

ThP-204 / DART IMAGING FOR MOLECULAR ACCUMULATIONS ON THE MUD SURFACE

Yayoi Hongo⁽¹⁾ - Motoshi Sakakura⁽²⁾ - Teruhisa Shiota⁽²⁾ - Kazumi Yoshiya⁽¹⁾

Tokyo Institute of Technology, Earth-Life Science Institute, Tokyo (1) - AMR, Inc., Mass Spec Application development, Tokyo (2)

Keywords: molecular mapping for non-structured surface, DART mass imaging, geochemical samples, on-surface derivatization

Introduction

Local accumulation of amino acids and oligopeptides on a surface of mud sediments was visualized using DART imaging. A movable stage and He cone were specialized to the analyses of unstructured samples. The image of surface distribution of the building blocks of biomolecules would help to see the correlation between local geochemical conditions and abiotic amino acid evolutions, suggesting the possible requirements of emergence of life on the primitive earth.

Methods

DART-OS (IonSense) was mounted on Q-TOF (Xevo QTOF, Waters) and Q-IMS-Q-TOF (Synapt G2, Waters). Activated helium supplied through a narrow hole (0.5 mm) ionized target molecules on a narrow surface, c.a. 2mm², of the non-structured mud samples collected at Ege-Jigoku in Beppu, Ohita, Japan. The sample stage was moving with a constant speed in 0.2 mm/sec. DART Source temperature was set at 430°C and derivatization was also examined using DBD-Py-NCS or TFATBDMS.

Results

Authentic oligoglycines, tetraglycine, glycilalanine, and alanilglycine dipeptides doped on the mud surface were detected with a spatial information. Heated Hegas dried and altered the wet mud surface locally, so that the position once ionized did not interfere the later ionizations at the next spot. The maximum sample size was limited to 3cm x 3cm due to the source configuration. The optimum distance between the sample surface and He cone was kept within 4mm. The position of the target molecule was reconstructed by the calculation using the constant moving speed of the stage and detection times of ions corresponding to the targets. On-surface chemical derivatizations with DBD-Py-NCS or TFATBDMS worked not only for enhancing the signals but also selecting the specific structures to the amino acids or oligopeptides. Gas-phase ion chromatography in the ion mobility mass spectrometry and MSn also enhanced molecular selectivity to characterize isomeric peptides having different amino acid sequence.

Conclusions

Unlike SIMS, MALDI, and DESI, DART ionization does not depend on the surface characters, such as electron conductivity or smoothness. DART ionization with the movable sample stage enabled us to visualize the accumulation of amino acids and oligopeptides on the unsmooth and heterogeneous surface of the unstructured mud sediment. On surface chemical derivatization was assisted by the ionization heat and profitable to detect the signals specific to the targets.

Novel Aspect

Amino acids and oligopeptide mapping of the heterogeneous surfaces of mud sediments was realized using DART imaging mass spectrometry.

ThP-205 / HIGH RESOLUTION TANDEM MASS SPECTROMETRIC ANALYSIS OF NATURAL PRODUCTS

Eun Sun Ji - Kwang Hoe Kim - Gun Wook Park - Jong Shin Yoo - Jin Young Kim

Korea Basic Science Institute, Biomedical Omics Group, Cheongju-si

Keywords: Natural product, High resolution tandem mass spectrometry

Introduction

Natural products are small molecules produced by organisms, including primary and secondary metabolites, which have complex chemical structures. Recently, advances in mass spectrometry (MS) with high sensitivity and specificity have enabled it to become a major tool in the analysis of natural products. However, the analysis of natural products is still challenge because of an insufficient reference data of tandem mass spectrometry. Therefore, to provide a sufficient reference is necessary for the characterization of various natural products.

Methods

We selected 1,400 natural products with the diversity of the structure. Then, we analyzed them using three different high resolution mass spectrometer including Orbitrap MS, Triple Time-of-flight (Triple TOF) MS and Quadrupole Time-of-flight (Q-TOF) MS. Tandem mass spectra were acquired at five different collision energies and they were used for the construction of spectral library.

Results

Currently, the library contains more than 20,000 tandem mass spectra for 1,400 compounds. Comparative analysis of tandem spectral data showed different fragmentation patterns according to the structure of natural products, instrument types and collision energies.

Conclusions

High resolution tandem mass spectral library of natural products was constructed and used to identify unknown compounds or secondary metabolites in complex sample.

Novel Aspect

A high resolution tandem mass spectral library for 1,400 compounds using authentic standards and purified compounds has been constructed.

ThP-206 / CHARACTERISING THE CATABOLISM OF PEPTIDES USING ION MOBILITY ENABLED HIGH RESOLUTION MASS SPECTROMETRY WITH MASS-METASITE INTEGRATION FOR DATA PROCESSING

Jayne Kirk⁽¹⁾ - *Ismeal Zamora*⁽²⁾ - *Antoni Riera*⁽³⁾ - *Tatiana Radchenko*⁽²⁾ - *Anna Escola*⁽³⁾ - *Russell Mortishire-Smith*⁽¹⁾ - *Yun Alelyunas*⁽⁴⁾ - *Mark Wrona*⁽⁴⁾

Waters Corporation, Wilmslow, UK (1) - Lead Molecular Design, S.L., Sant Cugat Del Valles, Spain (2) - Institute for Research in Biomedicine (IRB Barcelona), The Barcelona Institute of Science and Technology, Spain (3) - Waters Corporation, Milford, MA, USA (4)

Keywords: Mass-MetaSite, WebMetabase, IMS QToF, Biotherapeutics, ADME

Introduction

A growing need exists to characterize and optimize the absorption, distribution, metabolism and excretion (ADME) properties of biotherapeutics. The challenge is that few software packages available to characterise the clearance and metabolic fate of biotherapeutics. Here, we mine ion mobility high resolution mass spectrometric data for the analysis of biotherapeutic drug metabolism using the Mass-MetaSite and WebMetabase software platform for processing.

Methods

The studied peptides are 14-amino acid analogues of the hormone somatostatin. They were incubated in human serum at eleven time points between 0-48 hr. Data were collected on an ACQUITY UPLC coupled to an IMS QToF using HDMSE. Data were processed in a prototype version of Mass-MetaSite 5.1.9 able to read data directly from UNIFI 1.9.2 using the built-in Application Programming Interface (API) and also uploaded onto the server based application, WebMetabase 3.2.9.

Results

A number of catabolites for the parent peptides were identified over the time course and could be attributed to hydrolysis of the peptide backbone, with the main catabolites identified attributed to hydrolysis of the non cyclic portion of the peptides. A measure of the ratio of structurally matched to mismatched product ions found by Mass-Metasite provides confidence in catabolite assignment through data acquired by data independent acquisition (DIA). Approximately 99 % of the parent was rapidly (<10 min) turned over generating three key catabolites with minor metabolites, possibly attributed to hydrolysis of cyclic portion of the peptide, formed over the incubation time course. The use of HDMSE showed an improvement in data quality through alignment of precursor and fragment ions in both RT and DT (drift time) thus leading to improved spectral quality and an increase in identified fragment ions over the number of false positives for several peptides.

Conclusions

Mass-MetaSite and WebMetabase were used to characterize metabolic fate of the parent peptides showing that data can be collected routinely, processed and reviewed efficiently across time courses and treatments. A tightly integrated, UPLC HRMS informatics platform with high chromatographic performance, high mass accuracy and fast scanning capabilities afforded by the IMS QToF enabled rapid and confident detection of peptides and any corresponding catabolites.

Novel Aspect

Informatics and ion mobility HRMS workflow for detecting and reviewing peptide catabolites.

ThP-207 / SEMI-QUANTITATIVE LC/ESI/MS ANALYSIS USING PREDICTIVE MODELS OF ESI IONIZATION EFFICIENCIES

Jaanus Liigand⁽¹⁾ - *Piia Liigand*⁽¹⁾ - *Mari Ojakivi*⁽¹⁾ - *Karl Kaupmees*⁽¹⁾ - *Anneli Kruve*⁽²⁾

University of Tartu, Institute of Chemistry, Tartu (1) - Freie Universität Berlin, Institut für Chemie und Biochemie, Berlin (2)

Keywords: electrospray, ionization efficiency, prediction model, semi-quantitation

Introduction

As the omics are rocketing we gain insight into the chemical transformations occurring in living beings. The lack of authentic standard substances hinders obtaining deeper quantitative information about these processes. In recent years, artificial intelligence is improving different fields of life. We propose an approach of applying machine learning to overcome the lack of authentic standard substances via predicting the ionization efficiencies (IE).

Methods

Altogether, ionization efficiencies of 400 compounds were measured in electrospray ionization positive and negative mode in 29 mobile phases (altering the percentage of acetonitrile or methanol and covering a variety of the most used additives) in flow injection. For model development, descriptors from PaDEL and different machine learning approaches (multilinear regression, random forest regression and artificial neural networks) were used.

Results

For both ionization modes, universal models overall mobile phases were developed. For positive mode regularized random forest regression gave the best model with a residual standard error of 0.54 logIE units and R² 0.82. For negative mode, multilinear regression gave the best model with a residual standard error of 0.58 logIE units and R² 0.77.

The coefficients in these models are expected to depend on the instrument. Therefore, we propose using a set of calibrants to transfer ionization efficiency prediction models between instruments. This means analyzing a smaller set of calibration compounds (6) and samples together using the corresponding instrument and mobile phase in one batch. The relative IE values measured in this way can be used to fit a model between the logIE values observed for calibrants and descriptors. As a result, the fitting coefficient for the model are obtained which describe directly the system used to analyze the samples.

Validating this approach with gradient elution resulted in a 3.7-fold mismatch of concentration.

Conclusions

Universal ionization efficiency prediction models help to reduce the need for authentic standards with less than 4 times mismatch in concentration compared to a conventional comparison of peak areas that can result in a mismatch of several orders of magnitude. The approach only requires 2D structures of the compounds and a small set of compounds for system calibration.

Novel Aspect

Universal ionization efficiency prediction models are developed for positive and negative ESI, allowing semi-quantitative analysis on LC/ESI/MS with high accuracy.

ThP-208 / MSIML – AN OPEN-SOURCE PROGRAM FOR THE VISUALISATION AND PIXEL-BY-PIXEL CLASSIFICATION OF MASS SPECTROMETRY IMAGING DATASETS

Lukasz Migas⁽¹⁾ - Emrys Jones⁽²⁾ - Richard Chapman⁽²⁾ - Steven Pringle⁽²⁾ - Perdita Barran⁽¹⁾

Manchester Institute of Biotechnology, University of Manchester, Manchester (1) - Waters Corporation, Research Department, Wilmslow (2)

Keywords: Imaging Mass Spectrometry, machine learning, bioinformatics, DESI, software

Introduction

Mass spectrometry imaging (MSI) is used to interrogate the contents of biologically relevant samples often directly without sample preparation. The most readily employed methods (*i.e.* DESI and MALDI) generate vast amounts of data [1], often making the analysis cumbersome and time-consuming. In addition, the

resurgence of machine learning and deep-learning methodologies offer an opportunity to mine MSI datasets and provide a need for user-friendly and accessible software solutions.

Methods

MSIML was developed to enable rapid access to machine learning, clustering and multivariate methods for the analysis of MSI applications. The software suite is written in python and incorporates scikit-learn for machine-learning and neupy for neural-network training. Sections of whole mouse body, kidney and brain were analyzed using DESI on a Xevo G2-XS Q-ToF (Waters, UK) to demonstrate MSIMLs capabilities in the analysis of complex MSI datasets.

Results

Unsupervised multivariate analysis techniques is shown to greatly improve the interpretation of large MSI datasets, as well as aid differentiation of spectral and anatomical features. MSIML provides easy access to most commonly used methods including (*i.e.* PCA, NMF) as well as clustering (*i.e.* K-means, Birch) and embedding (*i.e.* Spectral Embedding, t-SNE or HSNE) algorithms. However the main strength of MSIML lies in its ability to interactively label MSI pixels of multiple files, enabling building of complex supervised classification models that can be easily applied to unknown samples[2]. Here, classification models were applied to segregate healthy and tumorous regions in mouse samples and differentiate tissue components. In addition, on-the-fly classification is demonstrated in a high-throughput environment where each mass spectrum is pre-processed and immediately labelled by previously build classification models, providing the user with more information and directions for further analysis.

Conclusions

Here, we present open-source platform (MSIML), an all-in-one program that can be used for the entire analysis workflow, including loading of single or multiple.imzML or Waters.raw files, pre-processing, spectral analysis, data mining, machine learning and sample classification. MSIML offers easy to operate GUI environment, enabling rapid interrogation of large datasets and providing valuable bioinformatics information.

Novel Aspect

Combined machine learning classifications and multivariate analysis of known and unknown imaging mass spectrometric datasets can enhance differentiation of complex samples.

References

1. McDonnell L., van Remoortere A, de Velde N., van Zeijl R., Deelder A., J. Am. Soc. Mass Spectrom., 21, 1969-1978 (2010)
2. Veselkov K.,Mirnezami R., Strittmatter N., Goldin R., Kinross J., Speller A., Abramov T., Jones E., Darzi A., Holmes E., Nicholson J, Takats Z., PNAS, 111, 1216-1221 (2014)

ThP-209 / BAYESIAN APPROACH TO AUTOMATIC MASS SPECTRUM PEAK IDENTIFICATION IN ATOM PROBE TOMOGRAPHY

Alexander Mikhalychev⁽¹⁾ - **Nikita Lappo**⁽¹⁾ - **Alex Ulyanenko**⁽²⁾

Atomicus OOO, Minsk (1) - Atomicus GmbH, Karlsruhe (2)

Keywords: atom probe tomography, Bayesian approach, automatic peak identification, peak labeling.

Introduction

Atom probe tomography relies on the analysis of mass-spectrometry data [1, 2] and requires preliminary identification of the mass-spectrum peaks. Establishing correct correspondence between the peaks in the

spectrum and the particular atomic species is crucial for accurate interpretation of the measured data [3, 4]. Manual peak labeling is a time-consuming procedure and vulnerable to errors [4]. We propose a Bayesian approach suitable for robust automatic peak identification.

Methods

We apply the Bayesian approach, previously developed for search and match analysis of X-ray powder diffraction data [5], to the problem of peak identification in atom probe mass spectra. The optimal peak association can be found as the maximal likelihood estimate of the specimen model for the given experimental mass spectrum. The designed technique is applied to the measured time-of-flight mass-spectrometry data.

Results

The Bayesian approach provides a general framework for derivation of physically grounded methods for peak identification. The initial information about the investigated specimen (*e.g.* presence or absence of some elements, possible correlations between presence of different ions, etc.) is encoded in prior probabilities for the sample models. The expected inaccuracies of the measurement and peak search are characterized by the likelihood function. Then, based on the measured data, posterior probabilities are calculated for the specimen models. The optimal peak labeling is provided by the most likely model (having the maximal posterior probability). By varying the probability models, one can construct different peak identification algorithms, suitable for different types of measured input data.

For testing, a peak identification algorithm based on a simple model of mutually independent inaccuracies of positions and intensities of the peaks was designed. Applicability and robustness of the algorithm was demonstrated for a number of experimental datasets.

Conclusions

The proposed Bayesian approach is universal and suitable for construction of different peak identification algorithms by taking into account different aspects of the data acquisition process. It is capable of correct analysis of overlapping peaks, of making model- and probabilities-based decisions when the mass-to-charge ratio is the same for different ions, and of taking into account possible global systematic errors (*e.g.* miscalibration).

Novel Aspect

Instead of “guessing” an automatic analysis algorithm, suitable for superseding manual peak labeling, we demonstrate derivation of such an algorithm from first principles.

References

1. D. J. Larson et al., “Local Electrode Atom Probe Tomography”, (Springer, New York) (2013).
2. Kelly T. F., Larson D. J., *Annu. Rev. Mater. Res.*, 42, 1 (2012).
3. Haley D., Choi P., Raabe D., *Ultramicroscopy*, 159, 338 (2015).
4. Hudson D., Smith G. D. W., Gault B., *Ultramicroscopy*, 111, 480 (2011).
5. Mikhalychev A., Ulyanenkov A., *J. Appl. Cryst.*, 50, 776 (2017).

ThP-210 / DEMYSTIFYING MATCH FACTORS IN LIBRARY SEARCHING

Arun Moorthy⁽¹⁾ - ***Anthony Kearsley***⁽²⁾ - ***William Wallace***⁽¹⁾

National Institute of Standards and Technology, Mass Spectrometry Data Center, Gaithersburg (1) - National Institute of Standards and Technology, Applied and Computational Mathematics Division, Gaithersburg (2)

Keywords: Algorithms, Electron Ionization Mass Spectrometry, Hybrid Search, Library Searching

Introduction

Match Factor is a term ubiquitous in mass spectral library-searching. It is a numerical measure of similarity between spectra used to sort a hit-list of potential identifications to a query. But what happens when multiple unique hits have very similar match factors? How do you decide which is correct? Here we detail how match factors used in library searching are computed [1,2], and present a decomposition method to help interpret hit-lists when ambiguities exist.

Methods

An EI spectrum can be represented as a vector where indices are the integral-mass of fragment-ions and values are the corresponding abundances. The masses of all peaks considered in a match factor computation is a set. This set can be decomposed into subsets: (1) matching peaks, (2) query only peaks, (3) library only peaks, (4) zero peaks. The ratio of sums of peaks in each subset can be used to interpret hit-list ambiguities.

Results

Simple and Hybrid match factors are computed as a function of peak intensities from the full set of peaks. This function is referred to as the modified dot product and differs from a standard dot product of two vectors in that the input values are rooted and the function output is squared and then scaled by 999. The difference between simple and hybrid match factors is in defining the subset of matching peaks. In hybrid match factors, this subset includes peaks from both spectra that differ by the molecular mass difference between the two compounds in addition to peaks that share the same mass.

When library spectra give the same match factors to a query, the ratio of sums from each subset of peaks can be used to discriminate hits. For example, a match factor of 800 due to peak intensity variability reducing the sum of matching peaks but with no unmatched peaks is more likely to be the correct identification of the query than the same score generated from a pair with many unmatched peaks in either the query or library spectra. Several example scenarios will be discussed.

Conclusions

Mass spectral library search algorithms use computed match factors to sort potential identifications of a query spectrum. When two or more hits have similar match factors, simply choosing the highest match factor may be incorrect. This presentation details a method of further decomposing match factors to assist analysts in this situation. An R implementation of this method is available through correspondence.

Novel Aspect

The decomposition of peaks into subsets is a novel method of analyzing match factors and greatly improves the interpretation of ambiguous library search results.

References

1. Stein S. E., Scott D. R., Journal of the American Society for Mass Spectrometry, 5, 859-866 (1994).
2. Moorthy A. S., Wallace W. E., Kearsley A. J., Tchekhovskoi D. V., Stein S. E., Analytical Chemistry, 89, 13261-3268 (2017).

For information please contact: arun.moorthy@nist.gov

ThP-211 / TOWARDS AUTOMATED REPORTING AND REVIEW OF METABOLIC PATHWAYS FOR LARGE MULTIMODAL MSI STUDIES

Teresa Murta - Alex Dexter - Spencer A. Thomas - Rory T. Steven - Chelsea J. Nikula - Efsthios Elia - Adam Taylor - Tingting Fu - Ala Al-Afeef - Bin Yan - Robin Philip - Kenneth N. Robinson - Josephine Bunch

National Physical Laboratory, NiCE-MSI, London

Keywords: MSI, automated data characterisation, metabolism, cancer

Introduction

Enormous effort has been made to understand the metabolic hallmarks of cancer [1]. MSI techniques are promising tools to visualise the spatial distribution of the molecules that constitute the metabolic pathways altered in cancer. Identifying and characterising which molecules can be measured is crucial to study their interactions. An automated tool to search the HMDB [2], collect and organise varied information regarding these molecules was implemented.

Methods

Either the list of detected masses or those identified as drivers of a multivariate analysis or machine learning classification result can be fed to this tool. The input is compared to a list of adducts generated from all molecules present in the database. For each potential assignment, information such as name, metabolic pathways in which it is known to be involved, and tissues in which it has been found, for example, are collected and organised for further exploration.

Results

This tool outputs lists of potentially measured metabolites or potential drivers of segmentation or classification results, as well as estimates of the metabolic pathway coverage (*i.e.* percentage of pathway constituents that can be measured), for example. Its outputs can be used to: (1) compare data acquired from the same sample but with different MSI techniques or with the same MSI technique but at different sites; (2) adjust the acquisition protocol to either maximise the number of measured molecules or to focus on a particular set of molecules (*e.g.* : those expected to be crucial to understand the particularities of metabolic alterations happening in tumours and surrounding tissues); and (3) restrict data post-processing analyses to particular sets of molecules. This tool will be part of an automated data processing pipeline developed to cope with the large amount of data generated by a 5-year multisite, multimodal, and multiscale project focused on improving our understanding of metabolic alterations occurring in various types of cancer.

Conclusions

This automated MSI data characterisation tool can be useful not only for optimisation of MSI acquisition protocols, but also to understand the usefulness of MSI techniques to investigate biologically relevant questions, in particular, those leading to a better understanding of cancer-specific metabolic alterations.

Novel Aspect

Automated collection and organisation of metabolic information available in MSI data, highly relevant to study the metabolic fingerprints of cancer using MSI data.

References

1. Hanahan D., Weinberg R. A., Cell, 144, 646-674 (2011)
2. Wishart D. S, Feunang Y. D, Marcu A., Guo A. C., Liang K., et al., Nucleic Acids Res., 46, D608–D617 (2018)

ThP-212 / MODPLUS: TOWARDS SEARCHING FOR HUNDREDS OF POST-TRANSNATIONAL MODIFICATIONS IN SHOTGUN PROTEOMICS

Seungjin Na - Jihyung Kim - Eunok Paek

Hanyang University, Dept. of Computer Science, Seoul

Keywords: bioinformatics, post-translational modification; software; tandem mass spectrometry

Introduction

Computational approaches to PTM identification using tandem mass spectrometry are severely limited when dealing with many types of modifications [1]. MODplus is a practical, unrestrictive modification search tool that works well for complete proteome databases, hundreds of modifications from the Unimod [2], multiple modifications in a peptide, partially tryptic peptides and so forth.

Methods

Sequence tag approach significantly filters out useless peptides from a database [3], and improved spectral alignment based on aligned multiple tags rapidly identifies top N candidate peptides. The spectral alignment defines regions (called gaps) within a peptide that were not aligned with any tags but possibly contain modifications, and the mass differences defined over the gaps can be interpreted by enumerating modifications in Unimod or user-defined list.

Results

We assessed the performance of MODplus with Q-Exactive data from human proteome (HEK293 cell lysate, ProteomeXchange: PXD001468) in terms of its search time and the number of identifications with varying the number of modifications. In case of a standard database search tool (MS-GF⁺ [4]), the search time scaled linearly with the number of input modifications and the number of identifications started to diminish with a dozen modifications. In contrast, MODplus search with ~700 variable modifications required only ~3-fold time than that with three modifications where MODplus search time was similar to that of the standard database search tool. Even as the search space got bigger with the number of modifications, MODplus lost almost no existing identifications (or identifications from a smaller search space) and yielded more identifications. We also demonstrated the utility of MODplus with TMT-labeled, phosphorylation enriched data (Human Skeletal Muscle, ProteomeXchange: PXD001543), which is of great interest in many laboratories.

Conclusions

A tag-based approach could resolve many computational issues in an unrestrictive modification search. MODplus successfully addressed both scalability and reliability challenges in searching for hundreds of modifications from human proteome datasets such as HEK293 cells and TMT-labeled phosphorylation enrichment. MODplus can be easily integrated with any search strategy or workflow to play its role as a part of a data analysis pipeline.

Novel Aspect

MODplus dramatically increases the depth of post-translational modification discovery by mass spectrometry-based proteomics.

References

1. Na S., Paek E., Mass Spectrom. Rev., 34, 133-147 (2015).
2. Creasy D. M., Cottrell J. S., Proteomics, 4, 1534-1536 (2004).
3. Mann M., Wilm, M., Anal. Chem., 66, 4390-4399 (1994).
4. Kim S., Pevzner P. A., Nat. Commun., 5, 5277 (2014).

ThP-213 / DEVELOPMENT AND VALIDATION OF SEMI-AUTOMATIC MALDI-HDX SAMPLE PREPARATION AND DATA ANALYSIS TOOLS

Sam MacIntyre - Thomas Nebl

CSIRO, Biomolecular Interactions, Parkville

Keywords: MALDI-TOF, Hydrogen-Deuterium Exchange, Automation, Biomolecular Interactions

Introduction

Use of MALDI-TOF as high-throughput hydrogen-deuterium exchange MS platform

Hydrogen deuterium exchange followed by Mass Spectrometry (HDX-MS) is a biophysical tool capable of probing protein/ligand interactions, conformational changes, and protein dynamics. Despite an increased number of applications, the expansion of the technology has been slowed by its intrinsic technical and analytical complexity (e . g. sequential online pepsin digestion and HPLC separation at pH 2.5, 0°C). Although most HDX-MS studies have been conducted with electrospray ionization (ESI), matrix-assisted laser desorption/ionization-time of flight (MALDI-TOF) mass spectrometry turns out to be a convenient tool for this purpose. MALDI-TOF combines the advantages of high speed of analysis and excellent sensitivity and accuracy of mass measurements with the capability to analyse peptide maps in a single spectral acquisition. With recent advancements in automated sample preparation entire HDX-MS experiments can be executed within one hour [1].

Methods

Development of an R-package for the semi-automated analysis of large MALDI-HDX datasets

We have developed an R-package to help analyse large MALDI-HDX datasets. This software facilitates the assessment of the quality of spectral segments for accurate centroid mass determination of identified peptides. It allows the user to semi-automatically validate, visualize and compare the relative deuterium incorporation of peptide features. The resulting HDX dataset can also be exported a format suitable for subsequent statistical analysis or 3D protein structure mapping using the MEM-HDX package [2].

Results

Development of a MALDI-HDX workflow for the analysis of peptide binding interactions

Coupling MALDI with HDX has been challenging because of undesired back-exchange reactions during sample preparation, spotting and analysis. In this report, we systematically evaluated different methods utilizing MALDI coupled with a commercially available high-throughput nanolitre dispenser capable of highly accurate and precise spotting of 96 samples in less than five minutes. We establish a sample preparation protocol using sub-zero temperatures and non-aqueous matrix solution that minimize back exchange. We applied our newly developed MALDI-HDX data analysis package to demonstrate specific binding of M2 anti-Flag antibody to FLAG or 3xFLAG peptides in a mix of 9 peptide calibration standards.

Conclusions

MALDI-HDX provides accurate linear epitope mapping data

The biophysical characterization of a monoclonal antibody recognizing different linear epitopes establishes the ability of MALDI-HXMS to provide accurate peptide binding data. This approach may be extended to glycopeptide or macro-molecular interactions, provided that peptide mapping is conducted carefully in order to minimize back exchange and confirm sequence coverage.

Novel Aspect

New software package facilitates the systematic analysis of MALDI-HDX data

We present a new R-package to help analyse large MALDI-HDX datasets. We used this software to systematically analyse peptide back exchange rates and other technical challenges associated with the coupling of MALDI and HDX. This has allowed us to establish a protocol utilizing a nanoliter dispensing system for rapid, automated MALDI sample preparation. This step is an essential pre-requisite for future use of MALDI-HDX as a high-throughput platform to study biomolecular interactions.

References

1. Pirrone G. F., Wang, H., Canfield N., Chin A. S, Rhodes T. S. Makarov A. A., Analytical Chemistry, 89, 8351-8357 (2017).
2. Hourdel V., Volant S., O'Brien D. P., Chenal A., Chamot-Rooke J., Dillies M. A., Brier S., Bioinformatics, 32, 3413-3419 (2016).

Daniel Ortiz - Laure Menin - Luc Patiny

EPFL, ISIC/SSMI, Lausanne

Keywords: MS tool, high resolution, protein fragmentation, internal fragments, contaminants, m/z deconvolution, teaching

Introduction

Cloud based services are becoming of main importance for many applications like collaborative working, large databases analysis and IoT (Internet of Things). Nevertheless, few scientific implementations are currently available online. For the last 10 years, we are developing core functionalities to allow mass spectra processing directly in the browser [1]. Recently, we compiled various tools to solve usual tasks encountered in MS research laboratories, available on the website <http://ms.cheminfo.org>. Those tools are open-source (MIT licence), don't require any installation and are always up-to-date (just need to reload the web page). In addition to classical tools such as calculation of isotopic distributions and simulation of fragmentation ions, the tool assists the user in finding possible molecular formulas from an experimental high resolution monoisotopic mass, associated with the list of PubChem candidates. Among others, the tool allows to check for the presence of known contaminants or compounds listed in a Google Spreadsheet that can be dynamically created by the user. This tool implements a matching algorithm that compares experimental and theoretical spectra in order to generate interactive reports that can be easily shared between colleagues. These algorithms were successfully used to analyze intricate systems involving for example protein metal complexes [2] or complex ionic liquid cluster deconvolution. <http://ms.cheminfo.org> also contains interactive teaching oriented tools that can be used in MS related lectures.

Methods

The libraries are implemented in JavaScript and available on GitHub. The main library can be downloaded on <https://github.com/cheminfo-js/molecular-formula>. The code is thoughtfully tested using jest and each commit is automatically processed in order to check for any broken code. Currently, there are over 150 tests that cover more than 80% of the code. Code syntax is enforced by over 200 eslint rules, that yields to a code easy to read, maintain and release under a MIT license.

Results

Several tools available on <http://ms.cheminfo.org> were unified and recently updated. These applications are daily used in our mass spectrometry facility of ISIC (EPFL), for a wide range of applications and projects. Among others:

Core function to generate isotopic pattern distributions. It can be calculated from a molecular formula or a chemical structure (Drawn or imported), with the possibility to select isotope abundances and charges. Based on the customizable FWHM, the system generates the theoretical isotopic distribution that can be compared with experimental data and exported as high quality figures (SVG).

Application to determine possible molecular formulas for a specific monoisotopic mass acquired by high resolution mass spectrometry (HRMS). It is possible to specify the molecular formula range for a single monoisotopic value but also to specify groups of atoms that, when combined, generate multiple combinations that will be matched and compared with experimental data. The system allows as well to specify the ion types (polarity, number of charges, adducts, etc.).

Easy to track classical contaminants inherent to the sample preparation such as detergents (Tween, SDS, Triton. . .), polymers (PEG, PPG. . .), contaminants from containers (phtalates, plasticisers) and solvent interferences. The experimental mass spectrum can be rapidly screened against the updated ms.cheminfo.org database of contaminants. Easycont browses the results ordered by similarity (%) and gives information about each identified contaminant.

Apm2s calculates theoretical fragment ions (a, b, c, x, y, z and internal fragments) from a given

protein/peptide sequence with any user defined modification such as post translational modifications, chemical modifications, non-natural amino acids, ligands or metal ions bound to the protein, and matches each individual theoretical isotopic pattern to the experimental mass spectra to provide a list of matches with similarity score. Very recently [2], this exhaustive Apm2s matching tool provided an abundance of data which gave critical insights into cisplatin-ubiquitin binding.

MS/MS fragmentation is a virtual and interactive fragmentation tool. Bonds can be selected from any chemical structure with a mouse click and all fragments and combinations of them are generated and compared to the experimental data.

The new deconvolution tool consider an experimental mass spectrum as a linear combination of spectra simulated from expected molecular formulas and this is achieved using SVD (Single Value Decomposition). This algorithm can deal with molecular formulas that yield to overlapping spectra, e.g. M^+ and M_2^{++} . For example it has been used to shed light on the different clusters and their relative intensity in a complicated mixture of ionic liquids.

Conclusions

Contrary to beliefs, JavaScript proves to be suitable for creating fast and efficient scientific tools from any Web browser and indeed are competitive with most vendor's applications. Many different analytical tools were gathered on the ms.cheminfo.org website such as isotopic pattern distributions generations, molecular formula finder from monoisotopic mass and several advanced tools to compare the experimental data to the automatically predicted isotopic distributions of thousands species. The applications associated with this automated workflow are endless and can certainly be of extremely valuable help for any user with internet access having complex high-resolution mass spectra data to analyze.

Novel Aspect

Complex problems involving MS measurements can not necessarily always be solved using commercial solutions. In top of that, commercial softwares are often expensive, inflexible, and their installation and upgrade procedures are generally time-consuming. In order to use the versatile ms.cheminfo.org toolbox, the only mandatory application is a web browser already installed on a computer. The free and open source applications of ms.cheminfo.org are in constant development, quickly evolving with the researchers requirements.

References

1. Patiny L., Borel A., J. Chem. Inf. Model, 53, 1223-1228 (2013).
2. Lee R., Menin L., Patiny L., Ortiz D., Dyson P., Anal. Chem. 89, 11985 (2017).

ThP-215 / COMPARISON OF HIGH- AND LOW-RESOLUTION MS DATA FOR DIRECT TISSUE PROFILING: MOVING FROM A LABORATORY INTO THE CLINICAL PRACTICE

Igor Popov⁽¹⁾ - *Anatoly Sorokin*⁽¹⁾ - *Evgeny Zhvansky*⁽¹⁾ - *Vasily Eliferov*⁽¹⁾ - *Alexander Vorobyev*⁽¹⁾ - *Stanislav Pekov*⁽¹⁾ - *Vsevolod Shurkhay*⁽²⁾ - *Alexander Potapov*⁽²⁾ - *Evgeny Nikolaev*⁽³⁾

Moscow Institute of Physics and Technology, Dolgoprudnyj (1) - N. N. Burdenko Scientific Research Neurosurgery Institute, , Moscow (2) - Skolkovo institute of science and technology, Skolkovo (3)

Keywords: clinical application, neurosurgical support, profiling.

Introduction

High-resolution mass spectrometers are widely used in research facilities around the globe, but applying mass spectrometry to clinical use requires reproducing the results on the low-resolution instruments. But in practice high- and low-resolution spectra obtained from one sample are not so similar as one could expect. In our work we analyze feature selection approach to determine features that are common for the different resolution spectra.

Methods

All experiments were performed on Thermo LTQ Orbitrap XL instrument using the direct method of extraction of lipids from brain tumor tissues followed by ESI. Multiple switching between Orbitrap and LTQ analyzers allows obtaining high-resolution spectra as well low-resolution ones almost simultaneously from the one sample.

Results

Both types of spectra were processed through the identical data analysis pipeline for denoising, aggregation, normalization, peak picking and peak alignment. Two datasets were then aligned to each other and analyzed. The low-resolution spectra contain about one-third of the peaks, detectable in high-resolution spectra; however, all major peaks are found in both types of spectra. Analysis of peaks selected as features to distinguish the tumor tissue from boundary brain tissue demonstrate that vast majority of selected features were preserved in both types of spectra. Moreover, we have created a mapping scheme, which allows us to use a classifier trained to separate tumor samples from non-tumor brain samples on high-resolution spectra with data from low-resolution spectra with a minor decrease in accuracy.

Conclusions

It was demonstrated that tissue profiling and identification techniques developed on the high-resolution mass-spectrometers could be transferred to low-resolution instruments that are commonly used in routine application.

Novel Aspect

It is shown that low-resolution spectra preserve distinctive features of brain tumor samples and could be used for sample classification in clinical settings.

ThP-216 / INTEGRATION OF HIGH MASS ACCURACY GC-TOFMS RESULTS USING EI AND FI IONIZATION METHODS

Yoshihisa Ueda - Ayumi Kubo - Kenji Nagatomo - Masaaki Ubukata - Takaya Satoh

JEOL Ltd., MS business unit, Tokyo

Keywords: GC/MS, EI, FI, Time-of-Flight Mass Spectrometer

Introduction

Electron ionization (EI), which accompanied by extensive fragmentation, is the most popular method used in gas chromatography-time of flight mass spectrometry (GC-TOFMS). The compounds are identified by mass spectral database search. On the other hand, the field ionization (FI) is also available in GC-TOFMS, which can obtain the molecular information from accurate mass. We have developed the software to integrate the GC/EI and GC/FI data obtained from GC-TOFMS.

Methods

The commercially available vinyl acetate resin was used as a sample for GC/EI and GC/FI analysis. The GC-TOFMS, (AccuTOF™ GCx-plus, JEOL) combined with a pyrolyzer (PY-2020iD, Frontier lab) was used for data acquisition. The thermal decomposition products were injected with split mode (100:1) and separated using GC capillary column DB-5msUI. The EI/FI combination ion source was used for GC/EI and GC/FI measurement.

Results

The data integration of GC/EI and GC/FI was performed with following steps. (i) The GC peaks were detected from GC/EI and GC/FI data separately. (ii) The mass spectra were extracted for each GC peak. (iii)

Database search was performed for extracted GC/EI mass spectra using NIST mass spectral library. (iv) The mass spectra from GE/EI and GC/FI data were linked using retention time information. (v) The accurate mass analyses were performed by using mass of molecular ion shown in GC/FI spectra and the masses of fragmentation ions in GC/EI mass spectrum.

The one useful point of the software is the confirmation of database search results with GC/EI mass spectra by molecular information of GC/FI spectra. The other is to support analysis of unmatched compounds by database search. Both elemental compositions of molecular ion and fragment ions, which represent substructure, can be estimated by accurate mass from GC/FI and GC/EI mass spectra, respectively.

Conclusions

We have developed the software to integrate high mass resolution GC/EI and GC/FI data. The functions are tested using GC/EI and GC/FI data obtained by Py GC-TOFMS analyses for vinyl acetate resin. The software is useful to overview the sample chemical information in complicated mixture.

Novel Aspect

Development of integration software for GC/EI and GC/FI ionization data acquired by high mass resolution time-of-flight mass spectrometer.

ThP-217 / BEYOND THE ACQUISITION OF FTMS SPECTRA - DATA INTERPRETATION AND VISUALIZATION OF ACCURATE DATA FROM COMPLEX SAMPLES

Alessandro Vetere - Wolfgang Schrader

Max-Planck-Institut für Kohlenforschung, Departement for Mass Spectrometry, Mülheim an der Ruhr

Keywords: FTMS – data analysis – complex mixtures –visualization - software

Introduction

The analysis of complex mixtures from energy related mixtures has been pushed forward during the last years by the development of powerful mass spectrometric machinery that is not only capable of producing ultra-high resolving and accurate mass data, but that is also of high sensitivity. Still problematic, is simple and coherent data analysis that allows to gather the relevant bits of data from the enormous variety of numbers acquired during a typical measurement. [1]

Methods

Vendor software that enables the complete analysis and non-targeted assignment of complex FTMS spectra from crude oil or biofuels, especially in a time-resolved manner, is to date not available. Therefore a mixture of third party software and self-written tools needs to be applied. An essential part during data analysis is adequate grouping and visualization of the results to convey the necessary information.

Results

The first step in data analysis typically involves the assignment of elemental compositions to each peak within a given spectrum. This step can be aided by the use of Kendrick mass scale[2] and the fact that within crude oil related samples one typically finds broad homologous series of compounds where the repeating units are either CH₂ and/or H₂. A first step in data accessibility involves grouping into so called heteroatom classes which define the amount of heteroatoms per molecule. For visualization and reevaluation of the data, as well as for the comparison of different datasets, several tools have been developed that first transfer the data into heavily structured MS Excel workbooks. From here graphs with varying degrees of detail can be easily created. An additional tool enables a straightforward overview of heteroatom class distributions by plotting the findings into a hexagonal grid. [3]

Most challenging is the handling of time resolved data for which a lab-built database approach was devised.

Conclusions

Mass spectra of crude oil related samples reach an average peak density of 100-200 signals per mass unit over a range of $100 \leq m/z \leq 1200$. For an aim-leading interpretation of such spectra not only mass resolution and accuracy are of utmost importance, also data interpretation and representation play a vital role. A broad and flexible set of visualization tools helps understanding the composition and the chemistry of such complex samples as crude oils.

Novel Aspect

New software tools were developed to ease up the data analysis and visualization pathway of complex crude oil samples.

References

1. Vetere A., Schrader W., ChemistrySelect, 2, 849-853 (2017).
2. Kendrick E., Anal. Chem., 35, 2146-2154 (1963).
3. Zhurov K. O., Kozhinov A. N., Tsybin Y. O., Anal. Chem., 85, 5311-5315 (2013).

ThP-218 / AUTOMATIC IDENTIFICATION OF METAL-BOUND BIOMOLECULES USING SNAP-LC, IMAGING, AND 2D-MS

Christopher A. Wootton⁽¹⁾ - *Matthew Willetts*⁽²⁾ - *Pui Yiu Yuko Lam*⁽¹⁾ - *Maria van Agthoven*⁽¹⁾ - *Mark P. Barrow*⁽¹⁾ - *Peter J. Sadler*⁽¹⁾ - *Peter B. O'Connor*⁽¹⁾

University of Warwick, Chemistry, Coventry (1) - Bruker Daltonics, Bruker, Billerica (2)

Keywords: FT-ICR MS, isotope analysis, LCMS, cancer therapeutics, data analysis

Introduction

The majority of chemotherapy utilizes metal-based compounds to preferentially attack cancer cells. Proteomic investigations of compounds and effects on cellular components is widely sought after to understand the mechanisms of action for potential drugs, but difficult to interpret. (1,2) To accurately identify and assign these species we present a new method to reveal metallodrug-modified species within complex data sets such as LCMS runs, MS imaging, and 2DMS.(3)

Methods

The SNAP Algorithm (V2.0), coded inside Bruker Data Analysis software (V4.0, 4.1, 4.2, and 4.3) was used to peak pick acquired data sets. SNAP-LC and MSi was coded using the internal visual basic environment inside DA to automate peak picking and processing.

All spectra were acquired using a 12T Solarix FT-ICR MS. LC-MS and MS/MS analysis was acquired using an EASY nLC II system coupled to a micro-ESI source, and used a water-CAN separation gradient and 15cm C18 nLC column.

Results

SNAP-LC analyses each and every scan from large LCMS acquisitions and produces a "metal chromatogram" of the specified metal searched. 3 This method was successfully applied to tryptic digest samples containing Pt (2), Os (4) and Ir metallodrugs. SNAP-LC was then applied to a proteomic investigation of an osmium metallodrug administered to a human ovarian cancer cell line, successfully identifying four modified drug target species within ~10,000 biomolecules observed, automatically and accurately.

SNAP-MSi reveals not only the location of all metal containing species (in a similar fashion to LA-ICP MS) but also the speciation information of the analytes detected, such as free drug identity or metal modified biomolecule, providing detail on drug targets and transformations. SNAP analysis of 2D data sets was of particular interest due to their complex nature and diagonal dependence of isotopes with respect to

precursor and data point location. (5)

Conclusions

Accurate analysis of complex data sets is the key to understanding large biological systems, but is hindered by sample complexity and the ability to differentiate species of interest. SNAP-LC, MSi, and 2D analysis allows us to exploit the unusual isotopic patterns of metal chemotherapy agents and accurately identify drug targets in complex samples regardless of prior identity. These techniques are also applicable to such as halogenated species, pesticides, etc.

Novel Aspect

Exploitation of metal isotope patterns using SNAP analysis reveals cancer drug biomolecule targets and drug transformations in biological samples.

References

1. Banerjee, S. ; Soldevila-Barreda, J. J. ; Wolny, J. A. ; Wootton, C. A. ; Habtemariam, A. ; Romero-Canelon, I. *Chem. Sci.* 2018, 9, 3177–3185.
2. Wootton, C. A. ; Sanchez-Cano, C. ; Lopez-Clavijo, A. F. ; Shaili, E. ; Barrow, M. P. ; Sadler, P. J. ; O'Connor, P. B. *Chem. Sci.* 2018, 9, 2733–2739.
3. Wootton, C. A. ; Lam, Y. P. Y. ; Willetts, M. ; van Agthoven, M. A. ; Barrow, M. P. ; Sadler, P. J. ; O'Connor, P. B. *Analyst* 2017, 142, 2029–2037.
4. Wootton, C. A. ; Sanchez-Cano, C. ; Liu, H. -K. ; Barrow, M. P. ; Sadler, P. J. ; O'Connor, P. B. *Dalt. Trans.* 2015, 44, 3624–3632.
5. Van Agthoven, M. A. ; Wootton, C. A. ; Chiron, L. ; Coutouly, M. A. ; Soulby, A. ; Wei, J. ; Barrow, M. P. ; Delsuc, M. A. ; Rolando, C. ; O'Connor, P. B. *Anal. Chem.* 2016, 88, 4409–4417.

ThP-219 / BACTERIAL WHOLE CELL TYPING BY MASS SPECTRA PATTERN MATCHING WITH BOOTSTRAPPING ASSESSMENT

Yi Yang⁽¹⁾ - *Yu Lin*⁽²⁾ - *Zhuoxin Chen*⁽³⁾ - *Tianqi Gong*⁽³⁾ - *Pengyuan Yang*⁽³⁾ - *Hubert Girault*⁽⁴⁾ - *Baohong Liu*⁽³⁾ - *Liang Qiao*^(3,5)

Department of Chemistry, Fudan University, Shanghai, China (1) - Research School of Computer Science, College of Engineering and Computer Science, The Australian National University, Canberra, Australia (2) - Institutes of Biomedical Sciences, Fudan University, Shanghai, China (3) - Laboratoire d'Electrochimie Physique et Analytique, Ecole Polytechnique Fédérale de Lausanne, Sion, Switzerland (4) - Department of Chemistry, Shanghai Stomatological Hospital (5)

Keywords: MALDI MS, bacterial typing, spectral pattern matching, bootstrapping

Introduction

Rapid and reliable bacterial typing is an urgent requirement for public safety [1]. Matrix-assisted laser desorption/ionization mass spectrometry (MALDI-MS) is now widely used to analyze bacterial samples[2]. Identification of bacteria at species level can be realized by matching the mass spectra of samples against a library of mass spectra of known bacteria. Nevertheless, identification accuracy should be further improved [3].

Methods

We propose a new framework for more accurate MALDI-MS based bacterial typing. This framework includes two steps. First, spectral pattern matching is applied to figure out the candidates of identification results, where we introduced two new measures for spectral similarity. Second, the identification results are further assessed by a novel bootstrapping model [4] to provide more reliable characterization.

Results

The performance of the framework was evaluated with a general dataset containing the mass spectra of 1741 strains of bacteria and another challenging dataset containing 250 strains, including 40 strains in *Bacillus cereus* group that were previously claimed impossible to be resolved by MALDI-MS [5]. Compared to the commonly used cosine correlation [6], our new similarity measures showed better performance. With the help of the bootstrapping assessment, the sensitivity of bacterial identification at the species level increased from 41% to 88% at the same error rate (5%) on the data set containing 1741 strains. Our approach also enabled to resolve the *Bacillus cereus* group.

Conclusions

This framework could provide much more reliable bacterial identification at both the genus and species level with the bootstrapping assessment. The method could promote MALDI-MS based bacterial typing and may further be applied to many other spectroscopy based analytics, *e.g.*, the MS/MS spectral pattern matching in proteomics and metabolomics.

Novel Aspect

To the best of the authors' knowledge, our method is the first to provide a statistical assessment to MALDI-MS based bacterial typing to improve accuracy and reliability.

References

1. Fitch J. P., Raber E., Imbro D. R., Science, 302, 1350-1354 (2003).
2. Chen J. H. K., Ho P. L., Kwan G. S. W., et al, J. Clin. Microbiol., 51, 1733-1739 (2013).
3. Mather C. A., Rivera S. F., Butler-Wu S. M., J. Clin. Microbiol., 52, 130-138 (2014).
4. Efron B., Ann. Stat., 7, 1-26 (1979).
5. Clinical and Laboratory Standards Institute, M58 methods for the Identification of Cultured Microorganisms Using Matrix-Assisted Laser Desorption/Ionization Time-of-Flight Mass Spectrometry, 1st ed. (2017).
6. Wan K. X., Vidavsky I., Gross M. L., J. Am. Soc. Mass Spectrom., 13, 85-88 (2002).

ThP-220 / DECOY METHODS FOR FALSE DISCOVERY RATE ESTIMATION IN HIGH MASS ACCURACY PEPTIDE SPECTRAL LIBRARY SEARCHES

Zheng Zhang - Meghan Burke - Yuri Mirokhin - Dmitrii Tchekhovskoi - Sanford Markey - Stephen E. Stein

National Institute of Standards and Technology, Mass Spectrometry Data Center, Gaithersburg

Keywords: peptide mass spectral library, target-decoy approach, PeptideProphet algorithm, false discovery rate, false positive

Introduction

Spectral library searching (SLS) is an attractive alternative to sequence database searching (SDS) for peptide identification due to its speed, sensitivity, and ability to include any selected mass spectra. While decoy methods for SLS have been developed for low mass accuracy spectral library, it is not clear that they are optimal to high mass accuracy spectra. Therefore, we report the development and validation of methods for high mass accuracy decoy libraries.

Methods

Reverse and Random decoy libraries are found to be suitable for high mass accuracy SLS. The first constructs spectra by reversing and the second by randomly replacing all peptide sequences in the library while maintaining the library's amino acid distribution. In both cases, the m/z 's of fragment ions are shifted according to the decoy sequences generated. No alteration to fragment ion intensities, peptide length, or precursor charge state are made.

Results

The utility of Reverse and Random libraries for target-decoy SLS in estimating false positives and FDRs was demonstrated using spectra derived from a recently published synthetic human proteome project. 1 For data sets from two large-scale label-free and iTRAQ experiments, these decoy building methods yielded highly similar score thresholds and spectral identifications at 1% FDR. The results were also found to be equivalent to those of using the decoy-free PeptideProphet algorithm. The Reverse method yields slightly higher thresholds presumably due to its keeping homology matches. However, differences in these methods are rather small. Using these methods for an example FDR estimation, MSPepSearch, which is a NIST-developed search software, led to 18% more identifications at 1% FDR and 23% more at 0.1% FDR when compared with other widely-used SDS engines coupled to post-processing approaches such as Percolator. An application of these methods for FDR estimation for the recently reported 'hybrid' library search method is also made.

Conclusions

The methods reported here for determining FDR values for high mass accuracy spectra enable results of library searches to be directly integrated with those of more widely used sequence searching methods. The availability of a reliable means for estimating these FDR values should encourage the use of SLS, which offers advantages in speed and sensitivity. It is shown that SLS can yield more peptide identifications than SDS, especially at the highest levels of confidence.

Novel Aspect

Reverse and Random decoy methods for high mass accuracy SLS permit the merging of SLS and SDS results, increasing assignment of peptides, leading to deeper proteome coverage.

References

1. Zolg D. P., Wilhelm M., Schnatbaum K., Zerweck J., Knaute T., Delanghe B., Bailey D. J., Gessulat S., Ehrlich H., Weininger M., Yu P., Schlegl J., Kramer K. ; Schmidt T., Kusebauch U., Deutsch E. W., Aebersold R., Moritz R. L., Wenschuh H., Moehring T., Aiche S., Huhmer A., Reimer U., Kuster B., Nat. Methods, 14, 259-262 (2017).
2. Burke M. C., Mirokhin Y. A., Tchekhovskoi D. V., Markey S. P., Heidbrink Thompson J., Larkin C., Stein S. E., J Proteome Res., 16, 1924-1935 (2017).

ThP-221 / TISANE DENOISING ALGORITHM FOR LARGE DENSE SPECTRA

Lionel Chiron ⁽¹⁾ - Marc-André Delsuc ⁽²⁾

Casc4de, Data analysis, Strasbourg, France (1) - IGBMC, Université de Strasbourg, Strasbourg, France (2)

Keywords: Denoising, large datasets, Petroleomics

Introduction

In Fourier Spectroscopies noise can critically hamper data analysis. In particular, for FTICR datasets characterized by a large number of points and in some cases, by a large density of lines, removing the noise is a non trivial task. Some denoising methods were devised (Chiron et al. (2014)) (Fabrice Bray and Rolando (2017)) able to copewith the size problem but for a large number of lines.

Methods

We present there a fast, paralisable algorithm TiSane which is able to tackle the denoising problem for both a large datasets and a large density of lines. This is done by performing the denoising by part, dividing the spectrum in small chunks on which the Sane algorithm misapplied. The resulting denoised chunks are then reassembled to produce the whole denoised spectrum.

Results

Recently, the invention of the ultra high resolution - harmonized FT-ICR cells (Nikolaev et al. (2011)) led to FTICR datasets containing millions of points. Here we show the processing of a petroleomic dataset obtained with such a cell, with $R \geq 500k$ at $m/z=400$, 4194304 points, with an estimated total number of peaks of about 15000. For Sane, with a computer having a RAM memory of 30Gb, the rank for this dataset cannot exceed 500. It is about 30 times under the value needed for the theoretical correct denoising with this method. Here TiSane got around the memory limitation by running with the following parameters: 4000 intervals, an adaptive rank and 15 processing cores. The processing took 40 seconds. Finally we observe that TiSane's robust cleanup gives access to the small isotopic patterns.

Conclusions

The processing with TiSane is quite fast considering the size of the dataset. It benefits from a double acceleration in comparison with the Sane algorithm by using both a divide and conquer approach and parallelization. Apart large and dense line data the algorithm can be applied to small datasets with the advantage of its speed.

Novel Aspect

TiSane can be also useful for Proteomics where denoising is crucial for peptides identification with automatic search engine like MASCOT.

References

1. Chiron, Lionel, Maria A van Agthoven, Bruno Kieffer, Christian Rolando, Marc-André Delsuc. Proceedings of the National Academy of Sciences of the United States of America 111,4, 1385–90 (2014)
2. Fabrice Bray, Lionel Chiron, Julien Bouclon, Christian Rolando. Anal. Chem. 89, 18589–18593. doi:10.1021/acs.analchem.7b01850. (2017)
3. Nikolaev, Eugene N., Ivan A. Boldin, Roland Jertz, Gökhan Baykut. Journal of the American Society for Mass Spectrometry 22, 1125–33. (2011)

ThP-222 / VALIDATION DATA PROCESSING BY ADVERSE 2.0 FREE SOFTWARE PACKAGE: ASSESSMENT OF THE ANALYTICAL REQUIREMENTS OF A GC-MS METHOD FOR THE DETERMINATION OF 16 PAHS IN SMOKED FISH

Arianna Piersanti⁽¹⁾ - *Tamara Tavoroni*⁽¹⁾ - *Luca Scrucca*⁽²⁾ - *Roberta Galarini*⁽³⁾

Istituto Zooprofilattico Sperimentale dell'Umbria e delle Marche "Togo Rosati", Sezione di Ancona, Ancona, Italy (1) - University of Perugia, Department of Economics, Perugia, Italy (2) - Istituto Zooprofilattico Sperimentale dell'Umbria e delle Marche "Togo Rosati", Centre for Development and Validation of Methods, Perugia, Italy (3)

Keywords: method validation, data processing, PAHs, smoked fish

Introduction

Polycyclic aromatic hydrocarbons (PAHs) are a wide group of carcinogenic environmental contaminants. Official methods are not set for food monitoring plans, therefore in-house developed procedures need to be validated to assess their compliance with the strict performance characteristics required by Regulation 836/2011 (1).

Methods

A new developed GC-MS procedure for the determination of 16-EU-PAHs in smoked fish was validated (2). The study was performed at three spiking levels (0.9, 3 and 7.0 $\mu\text{g}/\text{kg}$) in four different days for a total of 41

experiments. The data were processed with the in-house-developed package ADVerSE 2.0 (Analytical methods Validation results Speedy Elaboration) for the R statistical environment (3).

Results

In the European Union maximum levels are fixed in food for both, the sum of four marker PAHs, (benzo(a)pyrene, benzo(a)anthracene, benzo(b)fluoranthene and chrysene) and for benzo(a)pyrene alone (4). Linearity, precision, recovery, LOD and LOQ of a new method for the determination of PAHs in smoked fish were estimated using ADVerSE, a free software package developed by analytical chemists and statisticians. Measurement uncertainty was also calculated applying a top-down approach. The results were comparable with those obtained with the traditional spreadsheets, confirming that the method was fit for purpose (1,2). However, using ADVerSE, the Validation Report is produced in few seconds, overcoming tedious calculations and possible mistakes, especially when a large number of analytes is included in the method scope.

Conclusions

Thanks to the use of automated statistical tools to handle validation data, analytical chemists engaged in the official control can go back to their original job focusing the efforts to improve analytical protocols instead of studying possible approaches to uncertainty estimation or statistics methodologies to demonstrate the performances of their methods.

Novel Aspect

ADVerSE package allows a true feedback cycle going back and forward from method optimization to validation.

References

1. Commission Regulation (EU) No 836/2011. Official Journal of the European Union, L215, 9-16 (2011)
2. Lestingi C. et al. "A fit-for-purpose method to monitor 16-EU-PAHs in food – results of five years official control in two Italian regions". Food Additives & Contaminants Part A, 34: 1140-1152 (2018).
3. R Core Team (2018). R: A language and environment for statistical computing. R Foundation for Statistical Computing, Vienna, Austria. URL <https://www.R-project.org/>
4. Commission Regulation (EU) No 835/2011. Official Journal of the European Union, L215, 4-8 (2011)

ThP-223 / UNAMBIGUOUS PRECURSOR MASS ASSIGNMENT ENABLES ACCURATE AND SENSITIVE PEPTIDE IDENTIFICATION OF DATA-INDEPENDENT ACQUISITION MASS SPECTROMETRY DATA

Dong-Gi Mun - Sang-Won Lee

Korea University, Korea University, Seoul, Republic of Korea

Keywords: Data-independent acquisition, multiplexed MS/MS spectra, precursor mass, multi-stage database search

Introduction

Data-independent acquisition (DIA) approach has emerged as a promising strategy for robust and large-scale quantitation with the development of advanced mass spectrometer and data analysis tool [1]. However, accurate and sensitive interpretation of highly complex multiplexed MS/MS spectra acquired through DIA experiment still remains a challenge.

Methods

Novel DIA data analysis strategy are presented utilizing unambiguously determined precursor mass. Accurate monoisotopic masses of multiple precursor ions corresponding to multiplexed MS/MS spectra

were determined through mPE-MMR (multiplexed post-experimental monoisotopic mass refinement) procedure which was previously developed and applied for interpretation of multiplexed MS/MS spectra from DDA data [1].

Results

MS/MS data with determined precursor mass through mPE-MMR were searched against spectral library using in-house built spectral library search tool (The first stage search) under precursor mass and retention time tolerance. MS/MS data without peptide identification were searched against customized protein database using MS-GF⁺ search engine (The second stage search) [2]. The remained uninterpreted MS/MS data after the first and second stage search process were subjected to MODplus search engine[3] specified for unrestrictive modification detection. Up to 60,000 peptides were identified from 3hr gradient DIA experiments. Comparing to conventional DIA data analysis approach which does not consider precursor mass during database search, our strategy utilizing precursor mass effectively increase accuracy and sensitivity in peptide identification.

Conclusions

Interpreting multiplexed MS/MS spectra of DIA data were interpreted with assigned precursor mass through mPE-MMR and multi-stage database search. It provided increased sensitivity and accuracy in peptide identification.

Novel Aspect

Previously developed DIA data analysis tools considered elution time and fragment ions are considered for targeted extraction or spectral library search. There are few methods providing intuitive accuracy determinant parameter of identified peptides such as mass measurement accuracy.

References

1. Bruderer, R. et al., Mol Cell Proteomics., 14, 1400 (2015)
2. Madar, I. H. et al., Anal Chem., 89, 1244 (2017)
3. Kim, S. and Pevzner, P. A., Nat Commun., 5, 5277 (2014)
4. Na, S., et al., Mol Cell Proteomics., 7, 2452 (2008)

For information please contact: cacagong@korea.ac.kr

ThP-224 / CYCLOBRANCH: AN OPEN TOOL FOR ANNOTATION OF ISOTOPIC FINE STRUCTURES OF NATURAL PRODUCTS

Jiri Novak - Anton Skriba - Jakub Zapal - Marek Kuzma - Vladimír Havlíček

Institute of Microbiology, Czech Academy of Sciences, Prague 4, Czech Republic

Keywords: dereplication, natural products, isotopic fine structure, metabolomics, software

Introduction

CycloBranch (<http://ms.biomed.cas.cz/cyclobranch/>) is an open tool for de novo sequencing of non-ribosomal peptides and for dereplication of natural products from MALDI-MS, LC-MS, MSI (mass spectrometry imaging) and MS/MS datasets [1, 2]. Recently, a fusion of MSI data (imzML) with histology images was also reported [3]. Here, we present how CycloBranch can be used for more reliable dereplication of natural products benefiting from isotopic fine structure [4].

Methods

The high-resolution mass spectra of standard bacterial siderophores [5] were collected on 12T Solarix FTICR mass spectrometer (Bruker Daltonics, Billerica, USA) equipped with ESI source in positive ion mode. The results were verified by nuclear magnetic resonance (NMR) spectroscopy by ^1H NMR, ^1H - ^{13}C HSQC, ^1H - ^{15}N HSQC, ^1H - ^{13}C HMBC, TOCSY and J-resolved experiments (AVANCE III 700 MHz, Bruker Biospin GmbH, Rheinstetten, Germany).

Results

The fine isotope patterns of pyoverdines D/E - whose structures differ by (de)amidation only - were recorded with 840,000 @ m/z 1334 resolution in a conventional mass spectrum. If the minimum relative intensity threshold was 0.5%, the tool annotated 11 out of 12 theoretical peaks of pyoverdine D (m/z 1334.6) and 8 out of 9 peaks of pyoverdine E (m/z 1333.6). The fragmentation mass spectrum was collected with 550,000 @ m/z 1334 resolution. The complete m/z scale including the sequence-specific b-ion series was generated by CycloBranchin silico. The tool correctly annotated 201 out of 221 experimental peaks exceeding 0.5% relative intensity. Sixteen other peaks with intensities below 2.4% were correctly unannotated artefacts, and 4 remaining peaks with intensities below 1.6% were unknown. We observed 19 peaks with relative intensities 0.5-7.5% (without isotopologs) specific to pyoverdine E and only one peak m/z 601.24 (3%) unique to pyoverdine D. The remaining peaks corresponded to fragment ions of both pyoverdines.

Conclusions

The fine isotope feature of our open-source software called CycloBranch was reported on a showcase example of bacterial siderophores pyoverdines D/E whose masses differed by 1 Da. The quantitative results were confirmed by NMR analyses.

Acknowledgement

Ministry of Education, Youth and Sports of the Czech Republic (NPU LO1509) and Czech Science Foundation (16-20229S).

Novel Aspect

A free software utilizing the benefits of isotopic fine structure for dereplication of natural products from MALDI-MS, LC-MS, MSI, and MS/MS datasets.

References

1. Novak J. et al., J. Am. Soc. Mass Spectrom., 26, 1780-1786 (2015).
2. Novak J. et al., BBA - Proteins Proteom., 1865, 768-775 (2017).
3. Luptakova D. et al., Sci. Rep., 7, 16523 (2017).
4. Prichystal J. et al., Anal. Chem., 88, 10338-10346 (2016).
5. Pluhacek T. et al., Mass Spectrom. Rev., 35, 35-47 (2016).

ThPS-S12 - Nutraceutics directions

Chairs: Roberto Romero Gonzalez, Roberta Galarini

ThP-225 / MICROWAVE-ASSISTED EXTRACTION OF OLIVE LEAVES TO OBTAIN PHENOLIC COMPOUNDS WITH MODULATING ACTIVITY OF AMPK

Beatriz Martin-Garcia⁽¹⁾ - *Sandra Pimentel-Moral*⁽¹⁾ - *Ana Maria Gomez-Caravaca*⁽¹⁾ - *David Arraez-Roman*⁽¹⁾ - *Arsenio Muñoz de la Peña*⁽²⁾ - *Antonio Segura-Carretero*⁽¹⁾ - *Alberto Fernández-Gutiérrez*⁽¹⁾

University of Granada, Analytical Chemistry, Granada (1) - University of Extremadura, Analytical Chemistry, Badajoz (2)

Keywords: Microwave-assisted extraction, phenolic compounds, olive leaves, HPLC-ESI-TOF-MS, AMPK

Introduction

Phenolic compounds present in olive leaves possess beneficial effects to prevent chronic diseases [1]–[6]. Concretely, some of them provide the capacity to modulate AMPk activation, so they are proposed as therapeutic target in diseases associated with inflammation (obesity and diabetes)[7]. Microwave-Assisted Extraction (MAE) has demonstrated a great potential for obtaining extracts with higher content in these bioactive compounds.

Methods

An experimental design Box-Behnken has been performed in which temperature (50-150°C), solvent composition (15-80 % EtOH) and extraction time (5-20 min) have been evaluated to optimize the recovery of phenolic compounds that modulate the AMPk activation using MAE. The resulting extracts were injected into the HPLC-ESI-TOF-MS and phenolic compounds were analyzed following the method of Talhaoui et al. [8].

Results

A total of 47 phenolic compounds were identified in olive leaves extracts. Among them, a total of 6 compounds, that have demonstrated the capacity to activate AMPk [7], were quantified (hydroxyoleuropein isomer A, rutin, luteolin rutinoside, luteolin glucoside isomer A, hydroxyoleuropein isomer B and verbascoside). Thus, the total concentration of these phenolic compounds ranged from 0.060 mg g⁻¹ to 1.418 mg g⁻¹ depending on the different MAE conditions. The Box-Behnken analysis at 100°C, 5 min and 15% ethanol were the conditions that showed the highest concentration of phenolic compounds related to AMPK activation.

Conclusions

MAE has demonstrated to be a useful technique to improve the recovery of phenolic compounds able to modulate the activity of AMPk from olive leaves. Thus, 15 experimental conditions were assayed in the MAE experimental design in which temperature, solvent composition and extraction time were evaluated. The results showed that 100°C, 5 min and 15% ethanol were the optimal conditions to extract this kind of phenolic compounds.

Novel Aspect

These extract conditions could be used to obtain extracts for in vivo analysis and evaluate their possible anti-inflammatory effect (related to diseases such as the obesity and diabetes).

References

1. Benavente-García O., Castillo J., Lorente J., Ortuño A., and Del Rio J. A., *Food Chem.*, 68, 457–462 (2000).
2. Kontogianni V. G. and Gerothanassis I. P., *Nat. Prod. Res.*, 26, 186–189 (2012).
3. Susalit E., Agus N., Effendi I., Tjandrawinata R. R., Nofiarny D., Perrinjaquet-mocchetti T., and Verbruggen M., *Phytomedicine*, 18, 251–258 (2011).
4. Kontogianni V. G., Charisiadis P., Margianni E., Lamari F. N., Gerothanassis I. P., and Tzakos A. G., *J. Med. Food*, 16, 817–822 (2013).
5. Omar S. H., *Saudi Pharm. J.*, 18, 111–121 (2010).
6. Khalatbary A. R. and Zarrinjoei G., *Iran. Red Crescent Med. J.*, 14, 229–234, (2012).
7. Jimenez-Sanchez C., Olivares-Vicente M., Rodríguez-Pérez C., Herranz-López M., Lozano-Sanchez J., Segura-Carretero A., Fernández-Gutiérrez A., Encinar J. A., and V. Micol, *PLoS One*, 12, 1–22 (2017).
8. Talhaoui N., Gómez-caravaca A. M., León L., De R., Segura-carretero A., and Fernández-Gutiérrez A., *LWT - Food Sci. Technol.*, 58, 28–34 (2014).

ThP-226 / DEVELOPMENT OF ISOTOPE DILUTION-LIQUID CHROMATOGRAPHY/TANDEM MASS SPECTROMETRY FOR THE ACCURATE DETERMINATION OF TRANS- AND CIS-VITAMIN K1 ISOMERS IN FOOD

Byungjoo Kim⁽¹⁾ - ***Hyeyoung Lee***^(1,2) - ***Hong Hee Lee***⁽¹⁾

Korea Research Institute of Standards and Science, Center for Analytical Chemistry, Daejeon (1) - Dong-Eui University, Division of Applied Bioengineering, Busan (2)

Keywords: trans- and cis- Vitamin K1, Food, Isotope dilution mass spectrometry, LC/MS

Introduction

Vitamin K1 exists mostly in the trans form in most of foods such as vegetables and unfortified dairy products. However, infant formula and cereals fortified with synthetic vitamin K1 contains both trans- and cis-isomers. To evaluate the nutritional status of general foods (fortified or not), trans- and cis-vitamin K1 must be measured individually because only the trans-isomer has biological activity [1]

Methods

Isotope-dilution mass spectrometry based on LC/MS (ID-LC/MS) for the accurate determination of trans- and cis-vitamin K1 was developed in this study. A C30 column was used to separate and individually quantify the trans- and cis-vitamin K1 isomers. Vitamin K1-d₇ was used as an internal standard. The matrix effect profiles were observed by post-column infusion experiment to optimize sample preparation methods and selecting ionization mode for LC/MS [2].

Results

Trans- and cis-vitamin K1 were well separated with the C30 column. The matrix effect profiles showed that APCI was less susceptible to matrix effects than ESI. Therefore, APCI was used to minimize the difference in the ionization efficiency of the target analytes and their isotope analogues. The validity of the method was evaluated by measuring gravimetrically fortified samples of infant formula, kimchi cabbage, and spinach. The measurement results (after subtraction of pre-existing levels) agreed with the fortified values within their uncertainties. The method showed good repeatability and reproducibility (less than 2 % of relative standard deviation of multiple measurements) for both trans- and cis-vitamin K1 in infant formula, with relative expanded uncertainties of approximately 5 %. For unfortified vegetables, its repeatability and reproducibility were less than 2 % for trans- vitamin K1, but around 6 % for cis-vitamin K1 as its abundance in unfortified food is much less than trans- vitamin K1.

Conclusions

An ID-LC/MS method was established as a candidate reference method for the accurate determination of trans- and cis-vitamin K1 in food. The validity test indicates that the method is of higher-order metrological quality in respect of its repeatability, reproducibility, and measurement uncertainty.

Novel Aspect

Compared to currently available methods, this study developed an ID-LC/MS method for the accurate determination of individual trans- and cis-vitamin K1 in food samples.

References

1. Lowenthal J., Vergel Rivera G. M., J Pharmacol. Experimental Therapeutics, 209, 330 (1979).
2. Lee S., Kim B., Kim J., J. Chromatography A, 1277, 35 (2013)

ThPS-S13 - Cultural heritage and archaeology

Chairs: Klaas Jan van den Berg, Giorgio Mellerio

ThP-227 / MASS SPECTROMETRY ANALYSIS OF MEDIEVAL ARCHAEOLOGICAL JEWELS

Monika Cechová - Lukáš Kučera - Petr Bednář

RCPTM, Palacký University, Olomouc

Keywords: mass spectrometry, archaeological findings, SERJs, medieval jewellery

Introduction

The history of jewellery is long and goes back many centuries. Slavic women had worn characteristic metallic S-shaped-end ring jewels (SERJs) and those jewels were usually found near their temple skull bones. Laser desorption/ionization mass spectrometry (LDI-MS) allows to effectively study organic compounds profile on the SERJs surface and helps us to understand the way they were used - whether as a part of a headdress or knitted directly into braids of hair.

Methods

Synapt G2-S high resolution tandem mass spectrometer (Waters,USA) equipped with vacuum MALDI ion source was used for all experiments. Organic matter of unknown origin adhered to SERJ surface (from Dětkovice archaeological locality, Moravia region, Czech Republic) was analyzed directly on SERJs surface. LDI-MS and mass spectrometry imaging experiments were done without application of matrix.

Results

MSI of fiber residues and fiber-like imprints present on SERJs surface revealed the saccharide chains in the form of sodium adducts. Strong signals at m/z 201.0024, 363.0647, 525.1129, 687.1681, 849.2052, 1011.2562 and 1173.2977 were found in positive ionization mode and ascribed to cellulose. Besides, signals at m/z 437.0921, 599.2493 and 761.1912 were observed with lower intensity. The latter can be explained by a cross link cleavage of glucopyranose ring [1] in sugar chain and further confirm the presence of polysaccharides. On the other hand, no signal of animal origin was observed on SERJs surface.

A distinct group of signals (*i.e.* m/z 242.0504, 445.1452, 648.2383, 851.3283, 1054.4178) corresponding well with *N*-acetyl-glucosamine chains was found in MS spectra obtained from a small piece of organic matter hidden in corroded surface of another SERJ. *N*-acetyl-glucosamine is a common part of chitin present also in beetle elytron. We hypothesize that those signals correspond with remains of a corpse-eating insect.

Conclusions

SERJs were attached to a headdress by a plant thread (morphology and MS data). *N*-acetylglucosamine signals revealed presence of insect body remains. Systematic analysis is now focused on characterization of fibers with animal origin (*e.g.* residues of animal skin).

Novel Aspect

Laser desorption ionization mass spectrometry allowed characterization of the organic residues deposited in corroded surface of medieval metallic jewels.

References

1. Li H., Bendiak B., Siems W. F., Gang D. R., Hill H. H, Analytical Chemistry, 85, 2760-2769 (2013).

Acknowledgement

Czech Science Foundation (17-17346S), Ministry of Education Youth and Sports, Czech Republic (CZ. 1.05/2.1.00/19.0377; LO1305), Palacký University Olomouc (IGA_PrF_2018_027).

ThP-228 / DETERMINATION OF THE WINE AND FOOD TRACES FROM ARCHAEOLOGICAL MATERIALS BY GC/MS METHODS

Iulia Farcas⁽¹⁾ - *Zaharie Moldovan*⁽²⁾ - *Nastasia Pop*⁽¹⁾

University of Agricultural Sciences and Veterinary Medicine, Str. Manastur 3-5, 400372 Cluj-Napoca, Romania (1) - National Institute of Research and Development for Isotopic and Molecular Technology, INCDTIM, Cluj-Napoca, Str. Donat 67-103, Romania (2)

Keywords: Archeology, Ceramics, organic residues, GC/MS, molecular markers,

Introduction

The organic residues of the ceramic objects contain biomarkers produced by biochemical transformations of their original content [1-4]. The characterization of the lipid fraction in archaeological material is of great importance to obtain information on uses and dietary of people who used this material. For obtaining precise results it is important to identify presence of molecular markers in ancient materials. The aim of this study is to analyze archaeological materials with sensitive and selective methodology involving also optimal extraction methods. The method of this study allows the identification of markers of food, wine and resins.

Methods

The identification of the organic residues in archaeological materials was performed following few stages: a) extraction of lipid; b) compounds derivatization and c) instrumental analysis. The lipids were extracted using chloroform– methanol (2:1) mixture and then the extracts were derivatised by *N,O*-bis(trimethylsilyl) trifluoroacetamide (BSTFA). The analysis of the samples was carried out using a GC/MS system. The mass spectrometer was operated in the electron ionisation mode (70 eV). The mass range was scanned in the range of m/z 50-650.

Results

A robust analytical method for determination of markers of wine, food and resins from archeological materials has been developed. The mass spectra of molecular markers are registered and discussed. The method proposed here will allow that the technique GC/MS to can be used extensively for the identification of lipids and wine markers, understanding the structure of the residues preserved in the samples.

Conclusions

Chemical characterization of the organic residues of ceramic vessels, contributes to a better understanding of their uses of the ancient people. Identification of chemical markers, in conjunction with a traditional study of the ceramics, plays a relevant role in establishing what the original content of a specific vessel was.

Novel Aspect

The method was used for characterization of ceramic samples from centre of Transylvania, near city Turda which was built by Romans under name Potaissa. The samples are originated from Roman period. This investigation is made for first time in Transylvania region. The analyses allowed us to differentiate between archaeological material devoted to the production or preservation of wine and those used to store the vegetable oil or animal food.

References

1. Pecci A. et al., Journal of Archaeological Science 40, 4491-4498 (2013)

2. Garnier N. and Valamoti S. M., *Journal of Archaeological Science* 74, 195-206 (2016)
3. Giorgi G. et al., *Journal of Archaeological Science* 37, 1453–1457 (2010)
4. Salvini L. et al., *J. Mass Spectrom.* 43, 108–115 (2008)

ThP-229 / STABLE ISOTOPE COMPOSITION OF THE MODERN AND FOSSIL FRESHWATER GASTROPODS: IMPLICATIONS FOR PALEOCLIMATE RECONSTRUCTIONS

***Antonella Macri*⁽¹⁾ - *Paola Iacumin*⁽¹⁾ - *Donatella Usai*⁽²⁾ - *Veerle Linseele*⁽³⁾ - *Sandro Salvatori*⁽²⁾**

Università di Parma, Dipartimento di Scienze Chimiche, della Vita e della Sostenibilità Ambientale, Parma (1) - Centro Studi Sudanese e Sub-Sahariani (CSSeS), Treviso (2) - KU Leuven, Laboratory of Biodiversity and Evolutionary Genomics, Center for Archeological Sciences, Leuven (3)

Keywords: Stable carbon and oxygen isotopes, aragonitic shell, Central Sudan, Holocene

Introduction

Rapid Climate Changes (RCC) that occurred during the Holocene are very important to understand climate changes currently underway (1). The characterization of these abrupt events allows the evaluation of ecosystems and ancient populations from a resilience point of view.

Methods

Here we investigate variations in $\delta^{18}\text{O}$ and $\delta^{13}\text{C}$ of Mesolithic archaeo-malacofauna (*Pila* sp.) found in Central Sudan (2) at high-resolution intervals along the growth axis of shell. Furthermore, in order to correctly interpret the data on the fossil specimens and validate the use of the molluscs in paleoclimate studies, modern molluscs were collected from two locations along the White Nile valley.

Results

Comparison of $\delta^{18}\text{O}$ and $\delta^{13}\text{C}$ in archaeological shell and in modern specimens revealed significant differences. The $\delta^{18}\text{O}$ archaeo ranges from -10.97‰ to -6.23‰ whereas $\delta^{18}\text{O}$ modern ranges from -1.17‰ to +5.13‰ suggesting climate change. Moreover, the lower $\delta^{13}\text{C}$ values for modern shells suggest a variation of carbon source in the aquatic system.

Conclusions

This study demonstrates that a drying process occurred over time along the Nilotic area. The low oxygen isotope values of Mesolithic samples indicate significantly less evaporation of the Nile water and/or a different rainfall amount and seasonal distribution in comparison to the current river. (3)

Novel Aspect

This study provides the first application of the use of *Pila* sp. as a proxy for paleo-climatic information.

References

1. Alley R. B., Agùstsdóttir A. M., *Quaternary Science Reviews*, 24,1123-1149 (2005)
2. Salvatori S., Usai D., Zerboni A., *The AfricanArchaeologicalReview*, 28,177-211 (2011)
3. Leng M. J., Lewis J. P., *Environmental Archaeology*, 21, 295-306 (2016)

ThP-230 / MASS SPECTROMETRY STUDIES OF THERIAC FROM “SPEZIERIA DI SANTA MARIA DELLA SCALA” (ROME): MODERN TECHNOLOGIES TO LEARN FROM THE PAST

***Luisa Mattoli*⁽¹⁾ - *Enrico Flamini*⁽¹⁾ - *Michela Burico*⁽¹⁾ - *Valentino Mercati*⁽¹⁾ - *Rita Pagiotti*⁽²⁾ - *Maria Luisa Vázquez de Ágredos Pascual*⁽³⁾ - *Pietro Traldi*⁽⁴⁾**

Aboca Spa Società Agricola, Area Ricerca Fitochimica, Sansepolcro (1) - Università degli Studi di Perugia, Dipartimento di Scienze Farmaceutiche, Perugia (2) - Universidad de Valencia, Department of Art History, Valencia (3) - Città della Speranza, Istituto di Ricerca Pediatrico, Padova (4)

Keywords: Theriac, composition, analysis, mass spectrometry

Introduction

Theriac, a thick sticky liquid medicine, until the mid 18th century was regarded as universal panacea and was used to treat infectious diseases or poisonings. Theriac contained up to seventy different often exotic ingredients. The production of a proper Theriac took months and it was supposed to be left to mature for years [1] [2]!

Here the first approach to study the composition of an ancient Theriac sample from the “Spezieria di Santa Maria della Scala” (Rome) is presented.

Methods

Different chromatographic methods principally based on mass spectrometry techniques have been used to achieve a comprehensive Theriac characterization. The organic polar compounds were analyzed by UHPLC - qTOF-MS, in high resolution (2GHz). The volatile organic compounds were analyzed by means of GC-TQ-MS, using Head Space sampling or direct injection after sample derivatization with BSTFA/Pyridine. The metal part of inorganic compounds was analyzed by means of ICP-MS, the anionic part by IC-ConD.

Results

Very sharp chromatographic fingerprints of the Theriac sample have been taken. After have studying molecular mass characteristic and the fragmentation pattern of the various molecular species, by matching data with a “in-house” specific library, different compounds have been identified. In several case high resolution accurate mass have been crucial. Obviously, the identification step has been very intriguing! Some marker compounds specific of medicinal plants have been correctly identified but many other are very difficult to recognize, as consequence of the transformations due to Theriac preparation as well as to the aging of the sample. Afterwards, using external standard regression curves, the identified compounds have been quantified. Glucose and fructose, typically should be from honey, are the most abundant compounds, followed by glycyrrhizin, and salicylic acid. A typical flavonoid liquiritin has also been detected, confirming the presence of liquorice. Calcium, potassium, sulfur and iron at concentration of 0.1-0.5% (w/w) have been found, while chloride has been found at a lower value.

Conclusions

More than thirty organic compounds and twenty elements have been identified, but further studies have to be done to better elucidate the composition of Theriac. The preliminary results showed how potent is the analytical approach used. By means of data obtained from a comprehensive characterization it is possible to suggest the presence of interestingly ingredients in the old Theriac sample. The rational study of Theriac is a challenge but can be convenient to suggest new solutions for the modern human beings health care.

Novel Aspect

For the first time the composition of Theriac has been investigated at a molecular level using different mass spectrometry techniques, including HR-MS.

References

1. Griffin J. P., British Journal of Clinical Pharmacology, 58, 317-325 (2004).
2. Karaberopoulos D., Karamanou M., Androutsos G., Lancet, 379, 1942-1943 (2012).

ThP-231 / MOLECULAR CHARACTERIZATION OF THE PREPARATIVE LAYERS OF XXIInd SARCHOPHAGUS, A MASS SPECTROMETRY APPROACH

Chiara Melchiorre⁽¹⁾ - *Georgia Ntasi*⁽¹⁾ - *Leila Birolo*⁽¹⁾ - *Paola Cennamo*⁽²⁾ - *Giorgio Trojsi*⁽²⁾ - *Giancarlo Fatigati*⁽²⁾ - *Angela Amoresano*⁽¹⁾ - *Andrea Carpentieri*⁽¹⁾

Department of Chemical Science, University of Naples Federico II, Naples, Italy (1) - Faculty of Letters, Università degli Studi Suor Orsola Benincasa, Napoli, Italy (2)

Keywords: cultural heritage, restoration, extraction, separation, mass spectrometry

Introduction

The molecular characterization of historical artifacts is becoming an increasingly important aspect for cultural heritage [1]. Nowadays all the developments in extraction, separation and analytical methodologies can be helpful for challenging tasks, such as the chemical reconstruction of the natural compounds used by artists [2]. Any of the informations about the chemical composition of binders, oils, varnishes, etc., play a fundamental role both for restoration purposes and for the reconstruction of historical informations.

Methods

The first step of our workflow is represented by an ammonia extraction of chemical compounds from two micro samples from the preparative layers and two microsamples from external layers of a wooden sarcophagus, which can be dated towards beginning of the XXIInd Dynasty.

The methodology used for the characterization of extracted molecules are GC-MS (for monosaccharides and lipids analysis) and high-resolution MS (for proteins identification).

Results

According to extraction procedure, we divided analytes in three categories: monosaccharides, lipids and proteins.

Monosaccharides analysis led to the identification of binders extracted from all the layers. We identified Galactose, Glucuronic Acid, Ramnose and Arabinose, which are the main constituents of Arabic gum.

We could identify palmitic, stearic and oleic acid. The relative ratio of identified fatty acid is for sure affected by aging; their co-presence though is compatible with the putative presence of egg [3]. This hypothesis is strengthened by proteins identification. We could in fact identify in the external layers (the ones connected with the decoration of the sarcophagus) Lysozyme, one of the most abundant protein in egg.

Altogether, our data suggest that the painting techniques adopted was tempera. Moreover, in the preparative layers, we could unambiguously identify Collagen; the presence of this protein can be related to the technique used for the building of the sarcophagus itself.

Conclusions

Numerous approaches are focused on single classes of compounds; some others [4] suggest the possibility of using a single extraction procedure followed by different analytical strategies. The rationale at the basis of our approach starts from the latter assumption. However, the interesting aspect of our approach is represented by the heterogeneity of unambiguously identified molecules, which led to the characterization of the different layers sarcophagus.

Novel Aspect

The ingenuity of our methodology relies in the use of powerful analytical methodologies (*i.e.* high resolution MS) in the workflow thus improving the quality of obtained results.

References

1. Vinciguerra R, De Chiaro A, Pucci P, Marino G, Birolo L. *MicrochemJ.*, 126, 341-348 (2016).
2. Romero-Noguera, Julio, Inés Martín-Sánchez, María del Mar López-Miras, José Miguel Ramos-López, Fernando Bolívar-Galiano. *Electronic Journal of Biotechnology* 13, 6-7 (2010).
3. Sunwoo, Hoon H., and NaiyanaGujral. *Handbook of Food Chemistry*. 331-363 (2015).
4. Bonaduce I., Cito M., Colombini M. P. *Journal of Chromatography A* 1216, 32, 5931-5939 (2009).

AUTHOR INDEX

A

- | | | | |
|------------------------------|----------------|-------------------------------|------------------------|
| Abad-Santos Francisco | 526; 550 | Aldini Giancarlo | 281; 1104 |
| Abbatiello Susan | 428; 1184 | Aldrich Colin | 187 |
| Abda Julia | 474 | Alechaga Elida | 554 |
| Abdelrahman Anas | 812 | Alelyunas Yun | 573; 1217 |
| Abdelrazig Salah | 743 | Alesso Enrico | 972 |
| Abe Hiroko | 891 | Alexander Cameron | 760 |
| Abe Yuichi | 510; 636 | Alexander Morgan | 496 |
| Abedi Ghazaleh | 355 | Ali Marina | 399; 400 |
| Abell Andrew | 121 | Alkhawaja Bayan | 539 |
| Abi-Ghanem Josephine | 716 | Alladio Eugenio | 372 |
| Abuawad Alaa | 759 | Allal Ahmed | 501 |
| Acar Cevik Ulviye | 822 | Allard Pierre-Marie | 175 |
| Achkasov Sergey | 736 | Allegrone Gianna | 1048 |
| Acquaro Vinicius Jr | 383 | Allen Kevin | 217 |
| Acquavia Maria Assunta | 732 | Allers Maria | 943; 1157 |
| Ács András | 1028 | Alley William | 574 |
| Adachi Jun | 263; 510; 636 | Allmaier Günter | 100; 691; 800 |
| Adamko Darryl | 217 | Allouche Abdul Rahman | 298 |
| Adamo Annalisa | 288 | Allouche Abdul-Rahman | 47 |
| Adamov Alexey | 1179 | Allsworth Max | 552 |
| Adams Chistopher | 69 | Al-Majidi Mohammed | 569 |
| Adhami Neema | 751 | Almeida Fernando | 437 |
| Adhikari Sarju | 126 | Almici Camillo | 253 |
| Afonso Carlos | 200 | Al-Natour Mohammad | 760 |
| Afshar Sepideh | 171 | Alon Tal | 310; 326; 1208 |
| Agazzi Silvia | 451 | Alonso David | 311 |
| Agostini Marco | 571 | Aloui Ines | 499 |
| Agtuca Beverly J. | 81 | Alpan R. Serdar | 511 |
| Ahadi Sara | 413 | Alpha-Bazin Beatrice | 654 |
| Ahkami Amirhossein | 658 | Alsufyani Sultan | 1018 |
| Ahmadi Shiva | 397 | Altelaar Maarten A. F. | 408 |
| Ahn Joong Kyong | 1027 | Altmann Friedrich | 1019; 1020; 1023; 1024 |
| Ai Wanpeng | 450 | Altomare Alessandra | 1104 |
| Airavaara Mikko | 444; 770 | Alving Anjali | 649; 1158 |
| Aizikov Konstantin | 224 | Alving Kim | 551 |
| Ajaero Chukwuemeka | 1058 | Alyu Tekes Fayza | 821 |
| Ajello Elena | 920 | Amalian Jean-Arthur | 674 |
| Akamatsu Shusuke | 1132 | Aman Ahmed | 490 |
| Akashi Satoko | 1033 | Amante Eleonora | 372 |
| Akashi Takahiro | 369 | Amaral Bruno Cesar | 775 |
| Akashi Takumi | 710 | Amaral Phellipe Honório | 826 |
| Al Bougha Salam | 812 | Amirav Aviv | 310; 326; 1208 |
| Al Iskandarani Mohamad | 531 | Amoresano Angela | 1032; 1244 |
| Al-Afeef Ala | 147; 924; 1222 | An Hyun Joo | 535; 659; 1027 |
| Alava Thomas | 350 | An Jung Eun | 1084 |
| Alazzo Ali | 760 | Anastasi Federica | 285 |
| Albast Noor | 812 | Andersen Lars | 309 |
| Alberici Rosana M. | 1036 | Anderson David | 917 |
| Alcaraz Christian | 950 | Anderson Tim | 973 |

Anderton Christopher R. 81
 Andre Jacques 463
 Andriolo Cyrus V. 1114
 Andriot Isabelle 33
 Andrisano Vincenza 562
 Angeli Sergio 362
 Angella Lucia 423
 Angeloni Simone 1096; 1206
 Anichina Janna 660; 661
 Anliker Sabine 1049
 Annunziata Maria Grazia 692
 Antal Borbála 470
 Antoni Shaula 37
 Antonini Lorenzo 946
 Anumol Tarun 604
 Anzillotti Luca 451
 Aoki Jun 625
 Applová Lenka 702
 Aprea Eugenio 136; 314; 370
 Arashida Naoko 693
 Araujo Flabio 189
 Arias Ramirez Daniel 761
 Ariaudo Alessandra 540
 Arkind Galina 780
 Armengaud Jean 654
 Armirotti Andrea 403
 Armitage Emily G 615
 Arnaudguilhem Carine 132
 Arnold Norbert 594
 Arraez-Roman David 1237
 Arrandale Victoria 616
 Array Tabiwang N. 1119
 Arrebola Rosalia 556
 Arrey Tabiwang N. 224; 410; 433
 Arrizabalaga-Larrañaga Ane 847
 Artavv Viatcheslav 311; 575; 593; 1072
 Asakawa Daiki 151; 809
 Ashcroft Alison 779
 Assaker Richard 15
 Assis Ana Cristina 662
 Astefanei Alina 226
 Astorga-Wells Juan 398
 Astrup Arne 250
 Atik Ahmet Emin 511; 952
 Audano Matteo 694; 734
 Auger Serge 889; 990
 Ault James 779
 Aureli Federica 212; 213
 Austin Daniel E. 114
 Auwärter Volker 663; 664
 Avantaggiato Giuseppina 1006
 Avellone Giuseppe 1204

Awad Hanan 217
 Ayala-Cabrera Juan Francisco 481; 576
 Azimi Ali 399; 400

B

Baba Takashi 806
 Bacanelli Gisele 189
 Bacchiocchi Simone 1108
 Bach Veronique 885
 Bache Nicolai 275; 286
 Bachman Martin 492
 Bączek Tomasz 380; 542; 720
 Bae Jung-Woo 372
 Baessmann Carsten 341; 854
 Baghdasaryan Astghik 454
 Bagheri Nejad Ramin 377
 Bagilet Florian 893; 894
 Bagnati Renzo 1076
 Bahadory Flair 695
 Bähre Heike 1170
 Bai Yu 25; 450; 965
 Baier Vanessa 562
 Bailey Aaron 776; 788
 Bailey Derek J 823; 1184
 Bailey Laura 102
 Bailly-Chouriberry Ludovic 893; 894; 895
 Bajoub Aadil 341; 854
 Baken Kirsten 1094
 Baker David 778
 Baker Laura 660
 Baker Paul 1120
 Bakulev Vasiliy 818
 Bal Wojciec 105
 Balato Giovanni 255
 Balboni Fiamma 255
 Balcerek Martyna 888
 Balcke Gerd 713
 Baldi Talini Alessandro 967
 Baldini Andrea 255
 Balog Julia 740; 875
 Baloglou Aristeidis 107
 Baltora Sylvie 885
 Bande Rieke 278
 Bandini Andrea 744
 Bang Geul 1061; 1146
 Banks Phine 745
 Bantscheff Marcus 41
 Baptistão Mariana 317
 Barath Peter 557
 Barbante Carlo 1093
 Barbaro Elena 1093
 Barbas Coral 737

Barberis Elettra	16; 418	Becker Thomas	820
Barbero Luca	540	Beckman Joseph	676; 1182
Barborini Emanuele	247	Bednář Petr	1240
Barbosa Sergio	139; 850	Bednařík Antonín	632
Barcaru Andrei	696	Bednarz Hanna	901
Barceló Damià	60; 1078	Bedont Stella	1099
Barden David	875	Beecher Chris	708
Bärenfänger Melissa	1031	Behrendt Rayk	242
Baril Patrick	179; 1047	Bekeschus Sander	636; 653
Barnes Alan	269	Belarte Franco Carme	230
Barnes Loïc	47	Belaz Katia Roberta	1111
Barola Carolina	995	Belenky Alexei	1138
Barran Perdita	67; 173; 208; 305; 331; 745; 793; 810; 1159; 1186; 1196; 1219	Belford Michael	428; 1184
Barrett David	496; 743; 759; 1152	Bell Matthew	102
Barrow Mark P.	75; 91; 183; 202; 222; 959; 964; 1230	Bellagambi Francesca G.	37
Barry Simon	147; 924	Bellina Bruno	173; 305; 810; 1159
Barry William E.	604	Belloni Walter	312; 698
Barsch Aiko	80; 549; 748; 767; 890; 1142	Belov Mikhail	156; 679; 781
Bartels Benjamin	626	Benavides Felipe	762
Barthélemie Morgane	1044	Benedetti Barbara	313
Bártl Josef	697	Béni Szabolcs	1038
Bartolini Erika	1034	Benicky Julius	1026
Bartolini Giuseppe	255	Bennadji Zhora	346
Bartolini Manuela	281; 343; 562	Ben-Nissan Gili	115; 777; 798
Bartolucci Gianluca	395; 840; 841; 846	Bensch Wolfgang	106
Bartosz Piotr	1025	Bentlage Herman	294
Barwa Erik	304	Bentley Liz	367
Barylyuk Konstantin	74	Bentley Mark	583; 602
Baschung Yannick	28; 536; 807	Bequet Frédéric	716
Bateman Nicholas W.	436	Berdeaux Olivier	1128
Batkulwar Kedar Balaji	603	Berden Giel	19; 294; 295; 725
Batov Ilya	987; 1010	Berdyshev Evgeny	1115
Battistel Dario	1093	Berezovski Maxim V.	807
Baumann Marc	247	Bergin Sarah	35
Baumeister Tim	907	Bergquist Jonas	523; 1077
Bautista Yesica	1043	Bergström Lind Sara	523
Bavaro Teodora	541	Bernabò Nicola	1098
Bayat Parisa	99	Bernardini Roberta	541
Bayçin Deniz	511	Bernhardt Oliver M.	250; 426; 431
Baygildiev Timur	819; 835	Bernier Laurent	635
Baykal Ahmet Tarık	414; 525	Bersenkowitsch Nina K.	20; 104
Beauchamp Jonathan	31	Berta Máté	340; 570
Beaumal Jean-Yves	834	Bertelkamp Cheryl	1050
Bech Thøgersen Janne	577	Bertolini Giacomo	247
Becher Dörte	267	Bertotti Elisa	540
Becher François	256; 440	Bertuzzi Michela	784
Becher Simon	194	Betta Emanuela	314
Beck Jonathan	602	Beusch Christian	398
Beck Scarlet	1188	Beveridge Rebecca	173
Beck Sebastian	942	Beyer Martin K.	20; 104; 107; 301; 304
		Biagini Denise	37
		Bianchi Federica	451

Bianchi Michele	512	Bölsker Stefan	632
Bianco Giuliana	704; 732	Bolzani Vanderlan da Silva	735
Biasioli Franco	34; 370; 873	Bomgarden Ryan	438; 792
Bichl Gerlinde	968	Bona Gianni	327
Bienz Stefan	176	Bonna Arkadiusz	105
Bier Daniele	189	Bonnamy Anthony	688
Bierstedt Andreas	484	Bonometti Elisabetta	863
Bieske Evan	303; 309	Bontempo Luana	235; 1015
Biffo Stefano	273	Borghini Francesca	869
Bigda Justyna	720	Borik Adam	453; 605; 1089
Bigler Laurent	176	Borovcová Lucie	491; 969; 1168
Bijzet Johan	17	Borton Chris	595; 1073; 1107
Bidl Wolfgang	401	Borysik Antoni	119
Biliaut Emeline	716	Boschmans Jasper	552
Bilkova Zuzana	425	Bose Utpal	389
Billing Anja M.	248	Boselli Emanuele	188
Bilsel Gokhan	774	Bosnjak Masa	527
Binici Burcu	774	Bossak-Ahmed Karolina	105
Binnarsley Cory	1018	Botitsi Helen	1109
Birkemeyer Claudia	337; 452	Botrè Francesco	142; 666; 877; 878; 881
Birolo Leila	227; 1244	Bottaro Christina	1057; 1083
Bischoff Rainer	265; 650; 794	Bottiroli Riccardo	314
Bishop Dezmond	567	Botzanowski Thomas	320
Bisutti Vittoria	975	Bouchaert Emmanuel	458
Black Connor	852	Bouclon Julien	966
Black Rachelle	305	Boudier David	1041
Blackstock Lindsay	1068	Boulanger Pascale	117
Blakney Greg T.	201	Boulanger Renaud	33
Blanch Joe	339	Boulanger Renaud	864
Blanksby Stephen	21; 125; 192; 687	Bourassa Sylvie	402
Blennow Kaj	249; 954	Bourgin Marc	601
Blin-simiand Nicole	1087	Bourgogne Emmanuel	880
Blockhuys Frank	937; 1116; 1118	Bouville Alyssa	1005
Blokland Marco	160	Bouwer Marc	29
Bloomfield Nic	614	Bouyssièr Brice	132
Bo Young Oh	1110	Bovo Giorgio	273; 923
Boccard Julien	889	Bowers Michael	567
Bocharov Konstantin	1209	Bowman David	616
Bodet Baptiste	977	Boyarkin Oleg V.	308; 309
Bodner Martina	362	Boyaval Fanny	923
Bodo Enrico	128	Boychenko Oleksandr	435
Boehler Marc	601	Boyd Gabrielle	581; 883
Boersema Paul	17	Boyer Sophie	895
Böhm Günter	768	Boyken Scott	778
Boije af Gennäs Gustav	814	Boyle Billy	552
Boito Marta	848	Bozzolino Cristina	372
Bojko Barbara	65; 522; 528; 698; 712	Brabcová Ivana	1126
Bolfi Bianca	327	Braber Joris	1063
Boll Dmitriy	156	Braccia Clarissa	403
Bolla Jani Reddy	783	Brachtvogel Louise	404
Bolotnik Timofey	736	Bracini Chiara	228
Bolsheva Nadezhda	928	Bräcker Julia	287

Braconi Laura 840
 Braillard Stephanie 716
 Brakhage Axel A. 415
 Brandi Jessica 16; 288; 763
 Bray Fabrice 966
 Bredy Richard 298
 Bregy Lukas 216
 Brem Gottfried 441
 Brendle Katrina 296
 Breuer Matthew 599
 Brhozovskiy Alexander 292
 Bridge Wallace 735
 Brill Claudio 849
 Brina Daniela 273
 Brindle Kevin 147; 924
 Brink Andreas 935
 Brinkmalm Ann 249
 Brinkmalm Gunnar 954
 Brnakova Zuzana 1026
 Brochot Celine 885
 Brockmann Eike Ulrich 896
 Brockmeyer Jens 287; 863; 868; 876
 Broer Wim 997
 Broggin Massimo 699
 Brøndsted Nielsen Steen 309
 Bronova Irina 1115
 Bronzi Benedetta 849
 Broos Ludo 915
 Broughton Howard 171
 Brown Christopher 1035
 Brown Jeffrey 305; 593
 Brown Kerene 789
 Brown Michael 1134
 Brüchert Stefan 799
 Brückner Lea 711
 Bruderer Roland 45; 250; 426
 Bruderer Tobias 216; 454
 Brümmer Ina 863; 868
 Brunelli Laura 699
 Brungs Corinna 897
 Bruni Renato 1007
 Brünken Sandra 299
 Brunner Andrea M. 1050
 Brunner Andreas-David 767; 1188
 Bruno Giuliana 636; 653
 Bua Silvia 841; 846
 Bubis Julia 280; 405
 Bücher Katharina 805
 Buchmann William 499
 Buckley Tom 144
 Bucolo Claudio 374
 Budamgunta Harshavardhan 251

Budnik Bogdan 69
 Bugrova Anna 289; 292; 651
 Bukhari Tallat 89
 Bukowski Nick 339; 875
 Bull James 303; 309
 Bunch Josephine 147; 496; 804; 924; 1222
 Bunesco Andrei 716
 Buning Christian 549
 Buratto Steven 567
 Burhenne Jürgen 518
 Burico Michela 571; 1037; 1099; 1242
 Burke Meghan 1232
 Burla Bo 1148
 Burlingame Alma 42
 Burman Joachim 1211
 Burner Jake 690
 Burnum-Johnson Kristin 149
 Burobin Mikhail 324
 Burykina Anna 1174
 Busani Luca 1053
 Busch Florian 778
 Busch Michael 473; 581
 Buthier Pierre-Etienne 350
 Butt Craig 1073
 Buzy Armelle 556
 Buzzi Arianna 16; 418; 763
 Buzzigoli Emma 1053
 Buzzo Mauro 513

C

Cabaret Stéphanie 1128
 Caboni Maria Fiorenza 986
 Cabral Elaine C. 1036
 Cabrera Gabriela M 315
 Cabrices Oscar G. 660; 661
 Cai Chengyuan 337
 Cai Yizhi 1186
 Caiola Elisa 699
 Cakić Nevenka 700
 Calabrese Valentina 200
 Caleo Matteo 423
 Calleri Enrica 541
 Calvano Cosima Damiana 316; 578
 Camel Valérie 976
 Cameron Simon J. 193; 474
 Camilo Nilton Soares 827
 Camin Federica 235; 862; 1015; 1017
 Campbell Andrew D. 147
 Campmajó Guillem 139; 850
 Campos Oliveira Luan Felipe 317
 Camuto Cristian 881
 Can Nafiz Öncü 592; 822; 825; 832

Can Özge	511	Castellanos Pablo	109
Can Özgür Devrim	821; 825	Castellino Stephen	925
Cañal Villanueva María Jesús	419	Castro Gamboa Ian	735
Canavesi Rossana	514	Castro Valeria	1003
Candish Esme	698	Cataldi Tommaso	316; 578
Caneparo Valeria	16	Caterino Marianna	406
Cannizzaro Renato	703	Catracchia Marco	443
Cano Leticia	171	Catteuw Amelie	1006
Cao Jie	579; 677	Cauchi Michael	367; 1210
Cao Zhengwen	205	Caulfield Michael	285
Capitani Donatella	1105	Cautereels Julie	937; 1116; 1118
Capitoli Giulia	271; 916	Cavaglioni Alessandro	869
Capolongo Francesca	975	Cavalca Viviana	753
Capone Simonetta	366	Cavallarin Laura	1100
Capozzi Maria	578	Cavallo Gianni	674
Capozzi Martina	877	Cavanna Daniele	851
Cappellin Luca	515; 580	Cavasso Filho Reinaldo	307
Cappiello Achille	457	Cavazzana Annachiara	370
Caprioli Giovanni	24; 1096; 1102; 1206	Cazals Guillaume	864
Caprioli Richard	8; 122	Cazenave-Gassiot Amaury	1148
Capriotti Luigi	1034	Cazzamalli Samuele	520
Capuccella Marinella	1056	Ceccacci Sara	532
Capuccini Marco	1211	Cecchi Rossana	451
Caputo Anna Maria	669; 670	Cecchini Marco	285; 423
Carayol Jérôme	250	Cecconi Daniela	16; 288; 763
Carazzone Chiara	762	Cechova Eliska	38
Carboni Chiara	133; 824	Cechová Monika	1240
Cardasis Helene	701; 823	Čečková Martina	544
Careri Maria	451	Cederquist Henrik	113
Caretta Andrea	364; 370; 873	Celia Christian	531
Carita Correra Thiago	1059	Cemazar Maja	527
Carli Fabrizia	1053; 1117	Cennamo Paola	1244
Carlin Silvia	705	Cerli Chiara	97
Carlos Mickaël	678	Černigojc Urh	343
Carmany Daniel	473; 581; 883	Cerny Jiri	10
Carnevale Neto Fausto	735; 1160	Cerrato Giuseppina	863
Carnielli Virgilio Paolo	39; 520; 744	Cervi Gustavo	300; 302
Carollo Stefania	744	Cervinkova Barbora	391
Carone Matteo	1211	Cevenini Armando	406
Carpentieri Andrea	1032; 1244	Cha Hoon Suk	1027
Carradori Simone	531	Chacón-Patiño Martha L.	78; 201
Carrasco-Pancorbo Alegría	341; 854	Chae Won-Ki	560
Carrascosa Eduardo	303; 309	Chagovets Vitaliy	459
Carro Stephen	742	Chakraborty Asish	517
Caruso Donatella	219; 694; 734; 753	Chambert Stéphane	47
Caruso Michele	548	Chamoun Sherley	430
Carvalho Paulo Costa	775	Champy Pierre	1144
Casalone Cristina	1104	Chan Tak Wah Dominic	85; 1185
Casanal Ana	781	Chan Wan	748
Cashman Kevin D.	186; 1097	Chandler Courtney	1123; 1135
Casiraghi Costanza	520	Chang Eun Jung	981
Castagner Roberto	1106	Chang Wei-Hung	638

Chapman Richard	1219	Choi Hyun-cheol	981
Chapman Terry	539	Choi Inseoung	1195
Chardon Karen	885	Choi Jae	433
Charles Laurence	58; 547; 674; 1161	Choi Jung Hoon	1146
Chatelain Eric	716	Choi Min Je	372
Chatterjee Niladri	852	Chou Pei-Hsin	1052
Chau Man-Chu	447	Choudhary Alka	238
Chen Chao-Jung	637	Chrastina Petr	697
Chen Chung-Yu	1113	Chrenková Lucia	702
Chen Emily	410; 435	Christidis Nikolaos	258
Chen Peng	959	Chu Ivan K.	129; 938
Chen Shu	1071	Chu Rosalie	149
Chen Weibin	517; 574	Chulovskaya Anna	882
Chen Yan	701	Chumakov Peter	280
Chen Yet-Ran	638	Chun Kyung-Soo	372
Chen Ying-Lan	638	Ci Chenggang	677
Chen Yinjuan	567	Cianfarani Stefano	1053
Chen Yuju	436; 1029	Cianférani Sarah	320
Chen Yu-Qing	1133	Ciantelli Massimiliano	37
Chen Zhinan	436	Ciccotelli Valentina	970; 972
Chen Zhuoxin	1231	Ciceri Elena	973; 974
Chen Zibo	778	Cicognola Claudia	954
Chendo Christophe	1161	Cillero-Pastor Berta	420
Cheng Sychyi	318	Cimino Nicola	619
Cheng Yu-Hsiang	646	Cinquepalmi Valeria	316
Chetwynd Andrew	209	Ciociaro Demetrio	1053
Chevallier Olivier	852	Ciriaci Martina	1003
Chevolleau Sylvie	1005	Ciric Andrija	853
Chiaberge Stefano	1147	Cirigliano Adriana	315
Chiarelli Laurent	418	Ciuti Elisa	133; 824; 831; 837
Chiasserini Davide	436	Civitareale Cinzia	613
Chiavarino Barbara	128; 244; 306; 810; 1105	Cizkova Dasa	15
Chikaoka Yoko	639; 643	Claassen Christin	135
Chim Yuen-Ting	474	Cladière Mathieu	976
Chimjarn Supansa	58	Claesson Carina	430
Chin Agnes	842	Clark Robert	619
Chinello Clizia	253; 271; 273	Clarke Nigel	285
Chiron Lionel	964; 1233	Clarkson Paul	556
Chirulli Vera	133; 824	Clarysse Christophe	516
Chiu Cookson K. C.	75; 91; 158; 959	Claude Emmanuelle	627
Chiuminatto Ugo	997	Clausen Per Axel	1062; 1193
Chiumiento Lucia	238	Clavier Hervé	940
Chmara Kamila	712	Clay Alex	844
Cho Eunji	319; 1051	Clayton Olivia	562
Cho Kyung Chul	1014	Clemmer David	195; 1030; 1035; 1175; 1178
Cho Seungil	714; 718	Clench Malcolm. R.	665; 929
Cho Wonryeon	290	Closs Brigitte	1041
Cho Y. E.	1014	Cochran Blake	221
Cho Yunju	319; 1051	Codling Garry	587
Choi Bo Ram	766; 1153	Cody Robert B.	236; 321; 455; 500; 1126
Choi Dojin	1153	Coelho Mirela B.	317
Choi Dong me	981; 1014	Cogo Paola	39; 520; 744

Cohen-Dvashi Hadas	798	Costanzo Michele	406
Cojocariu Cristian I.	973; 974	Costello Catherine	683
Colamonici Cristiana	878	Cottica Danilo	66
Colas Cyril	179; 1047	Courtemanche Gilles	716
Colby Sean	196	Cousins Lisa	467; 494
Cole Laura	665; 929	Coviello Donatella	704; 732
Cole Richard B.	99	Cox David M.	979; 1000
Coletti Cecilia	810	Cox Dezeræe	43
Coll Sergi	554	Cox Jürgen	427; 1188
Collard Marie	893; 894	Cox Roger D.	367
Collinge Emma	183	Crawford Elizabeth A.	456
Colomban Silvia	865	Creaser Colin	773
Colombini Maria Perla	228	Creskey Marybeth	407
Commandeur Jan	1159	Crestoni Maria Elisa	128; 244; 306; 810; 939; 1105
Compagnon Isabelle	165; 298; 1027; 1180	Creutzmacher Robert	805
Compagnon Isabelle	47	Criscuolo Angela	687; 1119
Compagnone Dario	582; 669; 671; 860; 1098	Crossett Ben	436
Comstock Kate	701; 823	Crotti Sara	571
Comte Rahel	1049	Croubels Siska	1006
Comunità Fabio	877	Cruciani Gabriele	925
Concina Bruno	298	Cubadda Francesco	212; 213
Conquista Michela	1003	Cubała Wiesław	542
Conrad Theresa	415	Cucchi Ulisse	548
Conrads Thomas P.	436	Cudjoe Erasmus	820
Constant Eric	109	Cummins Scott	389
Conte Eleonora	16	Cunsolo Vincenzo	657; 866; 1112
Conte Raffaele	1053	Curcio Christine	917
Coombes Steven	358	Curcio Davide	878
Coon Joshua	1184	Cuttano Armando	37
Coope David	900	Cuyckens Filip	63; 131; 584; 832
Cooper Ben	1018	Cvačka Josef	354
Coppola Lamberto	366	Cveticanin Jelena	115; 780
Cordeau Emmanuelle	132	Cypel Marcelo	65
Cordeiro Carlos	662; 955	Cyr Terry	407
Cordelle Sylvie	33	Cyrus Papan	708
Cordes Henrik	562	Czabany Tibor	968
Cordwell Stuart J.	436	Czaplewska Paulina	257; 417
Corilo Yuri E.	78; 201	Czech Hendryk	457
Corinti Davide	244; 305; 810	Czech Joerg	549
Cornil Jérôme	1164	Czentnar Zoltan	663; 664
Cornwell Owen	779	Czerwenka Christoph	162
Corona Aurélie	885		
Corona Giuseppe	364; 703	D	
Correani Alessio	520	Da Ros Alessio	705
Correig Xavier	751; 905	Da Silva David	179; 1047
Correra Thiago C.	300; 302	Da Silva Katyeny Manuela	124
Corthals Garry	226	da Silva Samantha	996
Costa Jean	361	Dabrowa Kajetan	572
Costa José Luiz	884	Daems Elise	12
Costa Rafael G.	771	Dahlberg Jeffery	658
Costa-Lotufo Leticia Veras	1160	Dairii Ekaterina	750
Costantino Robin	978		

Dalla Massara Lorenza	39	De Pauw Edwin	338
Dalla Pozza Elisa	763	de Petris Giulia	946; 947; 948; 951
Dall'Acqua Stefano	236	De Pinto Vito	657
Dall'Anese Riccardo	255	De Re Valli	703
Dall'Asta Chiara	198; 851; 1007	de Rijke Eva	97
Dalmia Avinash	820; 1055	de Rosso Mirko	1106
Damborsky Jiri	38	De Salve Ilse	540
Damen Mirjam	563	De Sangro Giovanni	669
Damiani Tito	1007	De Santis Michele A.	1112
Damoc Eugen	156; 224; 776; 781; 799	De Simone Angela	562
Dan Kisoon	275	de Souza Guillermo	345
Dando Ilaria	763	De Vita Pasquale	866
Dane A. John	321; 1126	de Voogt Pim	1094
Dang C.	429	De Winter Julien	1164
Danhorn Andreas	908	Deák György	59; 503
Danikiewicz Witold	572	Dearth Stephen	755
Dannhorn Andreas	898	Deb Sanjukta	963
Dantas Rodolfo	1111	Debrauwer Laurent	185; 978; 1005
Danza Giovanna	395	Debussche Laurent	556
Darcsi András	1038	Decarli Denise	1099
Darebna Petra	385	Decker Trevor	114
D'Aronco Sara	571	DeCoste Jared	473
Darrouzet Eric	179; 1047	Defoort Martial	117
Das Sunit	191; 490	Degani Genny	281
D'Ascanio Vito	1006	Degano Ilaria	228
Dasler Daniela	1200	Degliesposti Gianluca	781
Dasouki Majed	812	Dehghani Alireza	640
Datwani Sammy	322	Dei Cas Michele	374; 392; 739; 1138
Dauly Clairey	1119	Dei Silvia	840
Dávalos Juan	307	Deininger Lisa	665
Davidsson Pia	769	Del Gaudio Federica	408; 532
Davies Gareth	492	Del Grosso Ambra	423
Davoli Enrico	1076	Del Grosso Erika	512; 514
Dayon Loïc	250	del Olmo-Iruela Monsalud	1165
De Andrea Marco	16	DeLand Phillip	492; 556
De Angelis Elisabetta	865	Delaporte Grégoire	976
De Fabiani Emma	694; 734	Della Latta Veronica	1053; 1117
de Folter Stefan	489	Dellanna Frank	284
De Giorgio Fabio	881	Delsuc Marc-André	220; 964; 966; 1233
De Goede Stefan	600	Delvaux Cédric	338
de Graaf Erik	679	Demianova Zuzana .	413; 723; 724; 726; 741; 754
de Graaff Floris-Jan	904	Demond Paul	473; 581
de Haan Noortje	166	DeMontigny Pierre	551
De Jesus Inacio Luciana	975	Deng Hu	357; 359; 607
de Jong Felice A.	708	Denisov Eduard	156
de Koster Chris	97	Denti Vanna	253; 271; 916
de la Torre Xavier	142; 666; 877; 878; 881	Deriu Davide	1147
De Leoz M. Lorna A.	555	Dervilly-Pinel Gaud	1165; 1166
de Moraes Zamarion Vitor	1059	Dervishi Ina	291
de Oliveira Silva Ebenezer	1111	Desai Reena	375
De Palmas Anna	230	Deshpande Abhijeet Yashwantrao	603
De Paola Massimiliano	1076	Desjobert Jean-Marie	361

D'Esposito Rebecca	174	Dopfer Otto	306
Despre Victor	109	Doria Luisa	193
Deti Serena	669; 670	Dornan Gillian	781
Deuscher Zoé	33	Dorner Brigitte G.	440
Dexter Alex	147; 924; 1222	Dorneyei Agnes	1121
Dey Sudhansu K.	149	Dorrestein Pieter	582; 670
Dhummakupt Elizabeth	473; 581; 883	dos Santos Pedro Reck	65
Di Capua Angela	956	Dossin Eric	583
Di Carlo Claudia	763	Dossmann Héloïse	940
Di Carro Marina	313	Dostál Štěpán	491
Di Francesco Antonella	866; 1112	Dotz Viktoria	166
Di Francesco Fabio	37	Dowling Kirsten G.	186; 1097
Di Marco Pisciotano Ilaria	613	Doyle Rory M	393
Di Minno Alessandro	753	Drahos László	686; 1028
Di Ottavio Francesca	582; 670	Drechsler Nicole	1188
Di Palma Serena	706	Dreisewerd Klaus	87; 498; 632; 896
Di Stefano Vita	1204	Drewello Thomas	681; 945; 1200
Diamanti Irene	995	Dreyer Benjamin	612; 807; 867; 874
Diamond Miriam	616	Dreyer Matthias	549
Diana Eliano	863	Droit Arnaud	402
Diana Pierrick	583	Drum Chester L.	749
Diaz-Fernandez Paloma	474	Dryahina Kseniya	458
Didier P.	556	Du Ruijun	624
Diederich François	74	Dubbelman Anne-Charlotte	584
Dietrich Demian	456	Dubkov Mikhail	324
Dietrich Dörthe	914	Dubois Christelle	256
Dietz Sophie	707	Dudek Bettina	434
Dilillo Marialaura	679	Duewer David L.	555
Dillen Lieve	131; 584	Duffieux Francis	556
Dimovska Nilsson Kelly	933; 934	Duffy Sarah K.	186; 1097
Ding Caroline	823; 961	Dugo Paola	237
Ding Chuanfan	323; 567	Duncan Kyle	149
Diniz de Oliveira Karina	884	Dunning Caitlin	573
Diskin Ron	777; 798	Dunston Mark	1186
Dittrich Marcus	415	Dunyach Jean-Jacques	1184
Dmitriev Alexey	928	Duport Catherine	654
Dobias Radim	379	Đuriníková Erika	545
Dodge Jeffrey	171	Dutertre Quentin	602
Dodgen Tyren	29	Dzhatdoeva Ayshat	1010
Dohnal Ilse	968	Dzuman Zbynek	851
Dojahn Jörg	413; 521; 535; 543; 648; 708; 723; 724; 726; 785; 1120; 1170; 1173		
Doležal Karel	733; 932	E	
Doll Sophia	286	Easter Ashley	781
Döll Stefanie	707; 709	Easterling Michael	80
Dominguez-Medina Sergio	117	Eberlin Marcos N.	317; 689; 826; 827; 884; 1036; 1111; 1139
Domínguez-Vega Elena	529	Eckey Kevin	325
Donald William	154; 221	Edler Stefan	1091
Doneanu Catalin	517	Edvardsson Vidar O.	396
Dong Qian	89	Edwards Giles	1018
Donnarumma Danilo	1034	Efimov Alexey	1010
Donnarumma Fabrizio	86; 476	Egan Joseph M.	1160

Egea González Francisco Javier 979
 Egli Stefana 1083
 Ehkirch Anthony 320
 Ehlert Sven 457
 Ehrhardt Franziska 664
 Eiffert Lara 216
 Eijkel Gert 79; 629
 Eiriksson Finnur Freyr 396
 Ejsing Christer 79
 Eke Zsuzsanna 982
 Ekroth Susanne 610
 El Assad Dana 628; 899
 El Baba Tarick 195
 El-Aneed Anas 217
 El-Armouche Ali 242
 El-Baba Tarick 1035; 1175
 Elia Efstathios 147; 924; 1222
 Eliferov Vasily 40; 1227
 Eliuk Shannon 823
 Elliott Christopher 852
 Ellis Shane 79; 629; 631
 Ema Tadashi 448
 Emami Khoonsari Payam 1211
 Emsley Paul 781
 Engbrecht Kristin 658
 Engeser Marianne 152; 816
 Enghild Jan J. 654
 Engler Anna 216; 729
 Enjalbal Christine 132; 864
 Eno Nathan 619
 Eo Jae Yeong 332; 464; 466
 Eom Kwon Yong 597
 Eparvier Véronique 1044
 Ercolini Carlo 970
 Erdemgil Yiğit 511
 Eric Bonneil 44
 Ernberg Malin 258; 283
 Ernst Robert 1123; 1135
 Eros Nikolett 1121
 Erra-Balsells Rosa 564
 Esch Patrick 194
 Escobal Micaela 345
 Escola Anna 1217
 Esen Cemal 456
 Es-haghi Ali 365; 377; 1194
 Espada Alfonso 171
 Espinola Enrique 348
 Espinoza Ed 96
 Estdale Sian 67
 Etienne Chris 792
 Evangelista Sara 1108
 Evers Waltraud 559

F

Fabini Edoardo 281
 Fabris Daniele 174; 187
 Fabrizz Sodr  Fernando 124
 Faccioli L cia Helena 1137
 Faivre Val rie 256
 Fakler Bernd 401
 Falcetta Francesca 699
 Falck David 166
 Falkenby Lasse 275; 286
 Fang Ru 149
 Fang Sui 1085
 Fangmeyer Jens 1162
 Fanti Federico 671; 1098
 Far Johann 338
 Faraone Immacolata 238
 Farcas Iulia 1241
 Farewell Anne 933
 Farinola Gianluca 578
 Farmani Zahra 1191
 Farouk Jaber 531
 Farr  Marinella 60; 585; 1079
 Faserl Klaus 209; 647
 Fassy Florence 556
 Fatigati Giancarlo 1244
 Feasel Michael 883
 Fedorova Yana 292
 Feher Ioana 1074
 Feil Regina 692
 Feil Stefan 478; 1080
 Fel Anna 257; 417
 Feliciano Rodrigo 1189
 Fellers Ryan 956
 Fenaille Fran ois 256; 440
 Feng Bainian 829
 Fenner Kathrin 1049
 Ferguson Gavin 735
 Ferguson Michael 437
 Fernandes Anna Maria Alves de Piloto .. 826; 827
 Fernandez Eloy 238
 Fern ndez-Guti rrez Alberto 341; 854; 1237
 Fernandez-Lima Francisco 199
 Fernandez-Penas Pablo 399; 400
 Fernando Velasquez-Garcia Luis 157
 Ferramosca Alessandra 366
 Ferrari Giulia 161
 Ferraris Francesca 212; 213
 Ferreira Fernando 346
 Ferreir s Nerea 1124; 1143
 Ferrero Federica 327
 Ferro Gian Luca 984
 Ferry Isabelle 191; 490

Fewer David P.	765	Franc Vojtech	425
Fialkov Alexander	310; 326; 1208	Francesconi Kevin	614
Figl Rudolf	1019; 1024	Francese Simona	665
Filatova Daria	585; 1079	Francini Alessandra	1088
Filip Szymon	291	Francis Daniel	35
Filipe-Soares Renata	147	Francis Jo	844
Filippi Antonello	448	François Jade	1005
Fiore Marli Fátima	765	Frank Johannes	691
Fiorentino Francesco	783	Frankevich Vladimir	289; 459
Fiori Maurizio	613	Frankfater Cheryl	1150
Fiorino Giuseppina Maria	868	Franssen Maurice	127
Firth Ben W.	684	Franzke Joachim	811
Firth David	67	Franzoi Alessandro	1015
Fischer Markus	74	Fraschetti Caterina	448
Fitzgerald Christopher	813	Frazer Bill	604
Fjeldsted John	1182	Frederiksen Marie	1062
Flagella Zina	1112	Fresch Marion	868
Flamini Enrico	1099; 1242	Friesen Vanessa	1058
Flamini Riccardo	1106	Frölich Lutz	249
Flanagan Kieran	1018	Fryčák Petr	491
Fleck Stefanie C.	1143	Fu Hongzheng	959
Fleiner Julian	601	Fu Lijuan	1120
Fleishman Sarel	777; 780	Fu Pingqing	1167
Fletcher John S.	933; 934	Fu Tingting	147; 924; 1222
Flinders Bryn	629	Fuchigami Sotaro	1033
Flitsch Sabine	1186	Fuchser Jens	901
Florêncio Maria Helena	662; 1054; 1075	Fujii Makiko	493
Floris Federico	158	Fujii Masaaki	46
Fodaroni Giada	1099	Fujita Setsuko	1122
Foerster Kathrin Isabelle	518	Fujito Yuka	1122
Fogliano Vincenzo	34	Fujiwara Makoto	957
Folberth Julica	764	Fukusaki Eiichiro	727
Fomsgaard Inge	388	Fuller Daniel	1035
Fonslow Bryan	698	Fuoco Roger	37
Forcisi Sara	80	Fursova Anastasia	157
Foreman David J.	150	Furst Tomas	1213
Forester Craig	42	Furtak Jacek	712
Forleo Angiola	366		
Fornarini Simonetta	128; 244; 306; 810; 939; 1105	G	
Fornelli Luca	956	Gabelica Valérie	12; 200; 245
Forster Yvonne	176	Gachotte Daniel	1145
Fort Kyle	156; 781	Gadaj Anna	144
Foster Greg	137	Gaeta Carmine	704; 732
Foster Leigh	792	Gaetano Anna	1099
Fostner Shawn	117	Gaggini Melania	1053; 1117
Foti Salvatore	657; 866; 1112	Gagne Didier	132
Fouquet Thierry	500; 504; 506	Gaikwad Manasi	612
Fourasté Isabelle	1041	Gaisl Thomas	729
Fournier Isabelle	15; 148; 458	Gajadhar Aaron	410
Frache Gilles	628; 899	Gajek Magdalena	673
Fraioli Valter	670	Galarini Roberta	991; 995; 1056; 1211; 1234
		Galba Jaroslav	557

Galembeck Fernando	689	Gentry Matthew	900
Galewska Marta	542	Georges Béatrice	893; 894
Galimberti Stefania	271; 916	Gerace Enrico	372
Galleggiante Crisafulli Annamaria	972	Gerbaux Pascal	1164
Galli Ubaldina	514	Gerbig Stefanie	922
Gallien Sebastien	410; 435; 436	Gerdle Björn	258; 283
Gallo Pasquale	613	Gerssen Arjen	160; 1008
Galvin Karen	186	Gethings Lee	198; 719
Gambaro Andrea	1093	Geven Aysun	832
Gamble Heather	494; 1055	Geyer Philipp E.	286
Gamez Bobby	788	Ghader Masoud	365
Gamoh Keiji	710	Ghaemmaghami Amir	759; 760
Gampe Nóra	1038	Ghafouri Bijar	258; 283
Gan Zhiwei	337	Ghassempour Alireza ...	259; 355; 377; 384; 1199
Ganapathy Ramesh	438	Ghezellou Parviz	259; 384; 1199
Gandhi Tejas	250; 426; 431	Ghiasikhou Sahar	520
Gangarapu Satesh	127	Ghidoni Riccardo	374; 1138
Ganzella Marcelo	795	Ghimenti Silvia	37
Gao Liang	1148	Ghislain Lucien	322
Gao Yufeng	967	Ghosh Anirban	245
Garancini Mattia	271	Giacomello Pierluigi	946
Garavaglia Silvia	512	Giacomozzi Linda	309
García Alonso José Ignacio	419; 422; 737	Giambelluca Sonia	520
Garcia Patrice	893; 895	Giammarinaro Pietro Ivan	1034
García-Campaña Ana M.	1165; 1166	Giampà Marco	901
Garcia-Ruiz Eva	1186	Gianbelluca Sonia	39
Gareis Manfred	980	Giannetti Luigi	613
Gariglio Marisa	16	Giannossa Lorena Carla	1203
Garikapati Vannuruswamy	259	Gianoncelli Alessandra	784
Garimella Sandilya B.V.	196	Gibb Sebastian	630
Garrido Frenich Antonia	836; 855; 979; 1081	Gidden Jennifer	488
Garzoli Stefania	946; 947; 951	Giera Martin	904
Gasilova Natalia	159	Gieselmann Volkmar	421; 640
Gaspari Marco	443	Gili Marilena	984
Gasperi Flavia	136; 314; 370	Gilmore Ian S.	147; 924
Gastaldelli Amalia	1053; 1117	Gimbert Yves	940
Gatchell Michael	113	Ginsberg Howard	191; 490
Gates Paul	94; 1125	Giordana Alessia	863
Gatineau David	99; 940	Giordano Nadia	423
Gaugg Martin Thomas	216; 729	Giordano Silvia	699
Gauglitz Julia	582	Giorgi Gianluca	230; 869
Gavard Rémy	158; 222; 959; 964	Girard Maria Fernanda Cifuentes	768
Gavrilova Svetlana	292	Girault Hubert	1231
Gay David	147	Giri Sanjeev	603
Gehm Christian	181	Giuffrida Maria Gabriella	1100
Geib Timon	141	Giuliani Alexandre	499
Geisslinger Gerd	1124; 1143	Giuliani Marcella M.	1112
Geisslitz Sabrina	410	Giusepponi Danilo	995; 1056
Gembruch Oliver	901	Gizaw Elias	1029
Genangeli Michele	460	Gjelstad Astrid	845
Genazzani Armando	512	Glabasnia Anneke	583
Gentili Matteo	831	Glabonjat Ronald	614

Glaros Trevor	473; 581; 883	Greco Donato	1006
Gobom Johan	954	Greenhalgh C.	929
Godejohann Markus	204; 549	Gregori Adolfo	671
Godin Simon	501	Greig Michael	559
Godoy Cancho Belén	1095	Grevelding Christoph	922
Goedderz Markus	640	Griffoni Francesco	1003
Goedecke Niels	1188	Grignani Elena	66
Goetze Sandra	436	Grimmer Christoph	82
Gokcen Taner	774	Grinfeld Dmitry	156; 224; 590
Golbabanezhad-Azizi Ali-Asghar	1057; 1083	Gröger Thomas	82; 203
Goldman Radoslav	1026	Grolla Ambra	512
Goldman Scott	285	Gronbach Manuel	1039
Goletz Theresa	540	Grondin Christel	880
Goleva Elena	1115	Grones Peter	733
Gomez-Caravaca Ana Maria	1237	Groselj Ales	527
Gómez-Ríos Germán Augusto	65; 160; 383; 485; 486; 528	Grossert J. Stuart	684
Gong Tianqi	1231	Grosshans Pete	339
Gonzalez Isabel	171	Grotemeyer Jürgen	297; 479; 680; 1021
Gonzalez-Gago Adriana	422	Grüning Anja	1101
González-Menéndez Pedro	737	Grünwald-Gruber Clemens	1020; 1023
Goo Young Ah	291	Grutza Marie-Luise	107
Goodacre Royston	745	Gruz Jiri	1213
Goodlett David	1123; 1135	Grzetic Martens Josipa	439
Goodwin Richard J. A.	147; 898; 924	Gu Hawei	735
Goral Ilona	1025	Guarcini Laura	448
Gordiyenko Yuliya	781	Guarienti Michela	784
Goren Ahmet C.	774	Gucký Tomáš	544
Gori Marco	1091	Guéhi Gbogbri François	864
Gornik Tjasa	527	Guéraud Françoise	978
Gorshkov Mikhail	223; 280; 405	Guercia Elena	865
Gorshkov Vladimir	90; 280; 405; 884	Guérineau Vincent	1044
Gorst-Allman Peter Charles	344	Guerrera Ida Chiara	406; 408
Gorti Santosh Kapil Kumar	1120	Guillemot Jean-Claude	556
Goryńska Paulina Zofia	522; 712	Guillén-Alonso Héctor	489; 630
Goryński Krzysztof	65; 522; 528; 712	Gulde Rebekka	601
Gorzolka Karin	707; 713	Gülser İlgü Ece	511
Gosciny Séverine	163; 1171	Günal Mehmet Yalçın	414
Gosetti Fabio	327; 418	Gundlach-Graham Alexander	212
Goshawk Jeff	1212	Günther Detlef	212; 232
Goto Takayuki	1132	Gürel Büşra	414; 525
Gotti Clarisse	402	Gurke Robert	1143
Gottschalk Christoph	980	Gurov Victor	324
Götz Sebastian	1142	Gutknecht Piotr	720
Götze Michael	118; 787	Gutmann Rene	478
Gozzo Fabio Cesar	775	Guttenberger Nikolaus	614
Gozzo Sebastiano	395	Güven Gülipek	511
Grabic Roman	453; 605; 1089	Guy Philippe Alexandre	583
Grabicova Katerina	453; 1089	Guzman Monsalve Alexander	202
Grasso Marco	273; 923	Gwon Mi-Ri	714; 718
Gratteri Paola	841	Gygi Steven	428
Gray Nicola	615	Gyr Luzia	811

H

Ha Na	828	Harmon Taylor	102
Haapala Markus	814	Harper John	599
Haapaniemi Esa	765	Hart Claire	1134
Haas Albert	1150	Hartmann Laura	805
Haas-Baumann Naemi	454	Hartmann Rafael	523
Haddingh Marrit	265	Hartungen Eugen	478
Haddleton David	509	Harwood D. Tim	375
Haefeli Walter E.	518	Hase William L.	112
Haenel Lorena	260	Hasegawa Yoko	623
Hage Christoph	787	Haselberg Rob	57; 529
Hager Jörg	250	Hashimoto Masahiro	586; 1126
Hahnefeld Lisa	1124	Hashiro Shuhei	693
Hajslova Jana	598; 851; 983	Hattendorf Bodo	232
Hakala Jani	461	Hatters Danny	43
Hakkarainen Minna	451	Haufler Robert	614
Halgand Frédéric	1163	Hauke Frank	1200
Halim Mohammad	117	Haupt Alexander	401
Halin Emilie	1164	Haupt Caroline	795
Halvorsen Trine Grønhaug	378; 386	Hauschild Jan-Peter	687
Hama Takenori	342	Hausner Lucrezia	249
Hamar Matyas	698	Havelund Rasmus	147
Hambly Andrew	1125	Havlíček Vladimír	379; 491; 1168; 1236
Hamidi Hamideh	377	Hawkes Jeffrey A.	1077
Hamm Gregory	147	Hawkins Aaron R.	114
Hammel Yves-Alexis	977	Hayakawa Eisuke	715; 1167; 1214
Hamon S.	556	Hayashi Akio	1149; 1156
HAMPL Ales	38	Hayashi Masahiro	328
Han Dohyun	275	Hayen Heiko	1129
Han Guanghui	776; 788	Hayes Aoife	186; 1097
Han Jong Yoon	784	Hayward Nicholas K.	267
Han Ling	168	Haywood-Small Sarah	929
Han Myung Sub	1109	Hazenberg Bouke	17
Han Sang Yun	941	He ChaoHui	1215
Han Song-Hee	828	He Lidong	78
Han Sun Ho	1065	He Xiang	660
Handelsman David J.	375	He ZiXing	1215
Handler Anne Mette	903	Headley John	1058
Hands-Portman Ian	75	Heaton Cameron	665
Hanel Gernot	478; 1080	Hebert Alexander	1184
Hankemeier Thomas	584	Heck Albert J. R.	62; 156; 408
Hanot Vincent	1189	Hee Sun Jeong	1110
Hanrieder Jörg	84	Heeren Ron M. A.	63; 79; 83; 278; 460; 629; 631; 915
Hans Volkmar H.	901	Heidelberger Sibylle	521; 535; 543; 785; 1173
Hansbauer Eva-Maria	440	Heijs Bram	166; 904; 921; 923
Hansen Brett	114	Heikaus Laura	260; 411
Hansen Harald S.	906	Heiles Sven	194; 922
Happo Naohisa	957	Heiner-Fokkema Rebecca	696
Harapanahalli Praveen	627	Heinilä Lassi	765
Harat Marek	712	Heintz Dimitri	890
Harder Alexander	1119	Heintze Christoph	92
Harmange Guillaume	69		

Helle Niklas	297	Hofman Jakub	544
Hellenäs Karl-Erik	1012	Höglund K.	954
Hellwig Nils	77; 786	Hoke Joachim	272
Hemdrup Anne Louise	169	Holanda Pereira de Souza Diego	1059
Hendrickson Christopher L.	78; 201	Holčapek Michal	1126
Hendriks Lyndsey	212	Hollender Juliane	184; 601
Heninger Michel	208; 1087	Holliman Christopher L.	64
Hennermann Julia	536	Hollingsworth Joy	658
Hennig Kristin	716	Hollubarsch Tobias	922
Henschel Henning	814	Hollywood Katherine	1186
Hense Genia	241; 685	Holmes Paula	599
Henson William	1150	Hölscher Christoph	910
Hentz Sebastien	117	Hölscher Dirk	626
Heo Seok	597	Holst Stefanie	166
Heo Sung Woo	1109	Hong Seong-Kuk	372
Herbig Jens	478	Hong Yuning	43
Herburger Andreas	20; 104	Hongo Yayoi	1216
Herman Mélanie	163; 1171	Honscha Walther	452
Herman Stephanie	1211	Hopfgartner Gérard	93; 768
Hermans J.	265	Hopkinson Alan C.	129
Hermansson Martin	79	Horiyanto Elvi	394
Hermosilla Carlos	922	Horn David M.	269; 799
Hernandez-Alba Oscar	320	Horning Ole Bjeld	275
Hernández-Mesa Maykel	1165; 1166	Horovitz Amnon	780
Herniman Julie	335	Horvath Rita A.	695
Herodes Koit	351	Horvatovich Péter	265; 696
Herrera-Ubaldo Humberto	489	Hoshino Tadahiro	1156
Herting Katherine	438	Hou Guixue	436
Herve Marius	109; 298	Houel Stephane	788
Hesse Dörte	412	Howes M-J	1041
Hessler Gerhard	549	Hoyas Sébastien	1164
Heywood David	719	Hrbek Vojtech	598
Hidayah Siti Nurul	807	Hsu Fong-Fu	1150
Hieftje Gary M.	231	Hsu Hsu Chen	1021
Hieta Juha-Pekka	633	Hu Jundong	1069
Higo Daisuke	643	Hu Peifeng	1189
Hill Chris	781	Hua Lei	32; 52
Hill Russell	1135	Huang Guang	1068
Hinz Christine	615	Huang Shih-Pei	1021
Hirsch Andreas	945; 1200	Huang Taya	1063
Histel Karen	894	Huang Ting	17; 426
Hitzenberger Jakob F.	1200	Huba Anna K	475
Hixson Kim	658	Huebsch Matthew	413
Ho Jenny	435	Hughes Sarah	1058
Hobby Kirsten	615	Huguet Romain	269; 428; 956; 1184
Hodgkinson Jane	35	Huhmer Andreas	269; 433; 730; 731; 746; 755; 961; 1145
Hodík Jakub	697	Hui Renjie	829
Höehr Nelci Fenalti	826	Hull George L. J.	1097
Hoeng Julia	602	Hummel Thomas	370
Hoerning Ole	286	Hunter Christie	413
Hoffmann Jan	462	Hupfer Yvonne	434
Hoffrogge Raimund	404		

Huppertz Andrea	518
Huppertz Laura	663
Huq Mohammad	53
Hurkova Kamila	851
Husen Peter	79
Husser Christophe	935
Hutmacher Robert	658
Hutterli Manuel	580
Huybrechts Bart	985
Hwang Euijin	1109
Hwang Geum-Sook	1127
Hwang Hye-Jeong	1195
Hwang Taeik	597
Hyland Katherine C.	595

I

Iablokov Vadim	921
Iacobucci Claudio	118
Iacumin Paola	848; 1242
Iakab Stefania-Alexandra	905
Iannece Patrizia	704; 732
Iannella Loredana	878
Iannone Michele	666
Ibrahim Yehia M.	196
Ickert Stefanie	942
Ifa Demian	1103
Ihling Christian	787
Iida Miyako	494
Iinuma Yoshiteru	1167
Ilag Leopold	451
Ilgın Sinem	833
Ilina Irina	280
Ilitchev Alexandre	567
Ilker Un	774
Illarionov Boris	74
Illiano Anna	1032
Indelicato Sergio	393
Indeykina Maria	289; 292; 651
Inglese Paolo	898; 908
Inomata Satoshi	1167
Inoue Rintaro	802
Inoue Takahiro	1132
Inzalaco Giovanni	869
Ioannidis Kostas	963
Ionas Alin	587; 616
Iriti Marcello	739
Isaac Giorgis	1190
Ishii Kentaro	328
Isobe Toshiaki	9
Isoyama Junko	510; 636
Issart Ambre	501
Ito Takahiro	494

Iván Béla	59
Ivanov Daniil	1151
Ivanov Mark	223
Ivanov Vladimir	324
Iwamoto Shinichi	151; 949
Iwan Volker	1021
Iwase Hirotaro	891
Iwata Futoshi	54; 330

J

Jaanus Liigand	338
Jaber Malak	1152
Jacob Minnie	812
Jacobsen Nicklas R.	1193
Jacobsen Sophie Chakroun	906
Jae Myoung Oh	1110
Jaffuel Aurore	269
Jahn Olaf	412
Jakubowski Norbert	206
Jamin Emilien	978
Jamnik Polona	420
Jancalek Radim	38
Janfelt Christian	903; 906
Jang Kyungho	828
Jang Kyung-Soon	1061
Janicka Monika	830
Janssen Maurice H.M.	50
Jansson Christer	658
Jantunen Liisa	616
Januario Robin	307
Janulyte Aurika	463
Jardines Daniel	142
Jaroch Karol	522
Jasak Julia	413; 979
Jaubert Murielle	895
Jaumot Soler Joaquim	1136
Javahery Reza	467; 480; 494
Jeanne Dit Fouque Dany	686
Jeanne Dit Fouque Kevin	199
Jeanvoine Yannick	112
Jedrychowski Mark	428
Jemielity Jacek	245
Jensen Keld A.	1193
Jensen Kenneth	614
Jensen Poul Henning	169
Jeong Kahee	1016
Jeong Wanseop	332; 464; 466
Jervelund Mark	79
Jervis Adrian	1186
Jesse Stephen	497
Ji Eun Sun	1217
Ji Luo	521

Ji Ming-Chao 688
 Jia Mengxuan 778
 Jiang Ping 1068
 Jiao Xingchun 1071
 Jimenez Connie R. 436
 Jiménez Medina María Luisa 855; 979
 Jin Hwan Hong 1110
 Jin Wen 614
 Jirásko Robert 1126
 Joblin Christine 108; 299; 678; 688
 Jobst Karl 616
 Johansen Andreas Eske 442
 Johansson Malin 610
 Johnson Mark 18
 Jöhren Olaf 764
 Jokela Jouni 765
 Jolliffe Charles 480; 494; 1055
 Jones Christopher 1189
 Jones Emrys 627; 629; 900; 1134; 1186; 1219
 Jones Glen 509
 Jones Hugh 222
 Joos Thomas O. 13
 Joosten Marc 265
 Jordan Alfons 478; 1080
 Jordan Johannes 401
 Jordens Jan 57
 Jørgensen Thomas J. D. 169; 1202
 Jornet Rafael 230
 Josephs Jonathan 776; 788; 823
 Jung Eui Hyun 560; 561
 Jung John 1115
 Jung Wonhyeuk 155
 Jung Y. 290
 Jung Youngae 1127
 Junot Christophe 256; 440
 Jusko Pavol 299; 688

K

Kaabia Zied 893; 894
 Kaboord Barbara 438
 Kaczmarczyk Kinga 1025
 Kadesch Patrik 922
 Kaever Volkhard 1170
 Käfer Uwe 82; 203
 Kaftan Filip 907
 Kahen Kaveh 1060
 Kähler Kevin 1082
 Kailasam Srividya 870
 Kajita Ryo 639
 Kalenius Elina 208; 1196
 Kaling Moritz 596
 Kalli Anastasia 755

Kallinich Constanze 641
 Källsten Malin 523
 Kaltner Florian 980
 Kamgang Nzekoue Astride Franks 24; 1102
 Kamihoriuchi Bui 54; 330
 Kamiya Katsumasa 151
 Kammeijer Guinevere 166
 Kana Antonin 1070; 1198
 Kanawati Basem 80
 Kang Byeongwon 332; 464; 466
 Kang Dukjin 264; 642
 Kang Gihaeng 597
 Kang Guyoung 1051
 Kang Hee-Gyoo 262; 416; 717; 721
 Kang Hoil 597
 Kang Hye Soon 981
 Kang J. 290
 Kang Min-Jung 261; 534
 Kanicky Viktor 505
 Kaplancikli Zafer Asim 822; 833
 Kappes Manfred M. 296
 Karamessini Denise 58
 Karayel Başar Merve 414; 525
 Kargosha Kazem 365
 Karimi Kudratullah 462; 1040
 Karjalainen Maaret 669
 Karnoub Melodie Anne 15
 Karpov Dmitry 405
 Karras Gabriel 109; 298
 Karroca Ornela 1093
 Karst Uwe 325; 897; 914; 1162
 Kasahara Keiko 263
 Kashimura Akane 936
 Kashio Taiki 710
 Kaska Łukasz 720
 Katarzyna Pawlak 524
 Kato Daiki 1033
 Katz Aliza 777
 Kaufman Kimberley 399; 400
 Kaupmees Karl 1218
 Kauppila Tiina 633
 Kawamura Takeshi 639; 643
 Kawano Chihiro 715
 Kay Richard 391
 Kayser Sven 628
 Kazanc Emine 898; 908
 Kazemi Seyed Mahdi 259
 Kearsley Anthony 95; 1221
 Kee Chee Leong 842
 Kehrwald Natalie M. 1093
 Kéki Sándor 59; 470; 477; 503; 570
 Kelar Jakub 505

Kelleher Neil	956	Kim Jae Young	465; 910
Keller Sascha	272	Kim Jaehan	1027
Kellner Udo	901	Kim Jeong Ah	1146
Kelly Alan K.	186; 1097	Kim Jeongkwon	1022
Kelly Isabelle	402	Kim Jihyung	1223
Kelly John	716	Kim Jin Young	1217
Kelly Patrick	21	Kim Jin-Mun	988
Kempa Emily	331; 1186	Kim Jun Hyung	597
Kempf Jürgen	664	Kim Kwang Hoe	1217
Kempf Stefan	884	Kim Kwonseong	264; 642
Kempke Lisanne	19	Kim Kyoung Heon	1027
Kenwright Alan	358	Kim Min-Gul	828
Kepinski Stefan	733	Kim Seung Yong	332; 464; 466
Kernan Jeffrey	619	Kim Sook Heun	1109
Kero Frank	1091	Kim Suhkmann	766
Keshet Uri	310; 326	Kim Sunghwan	319; 1051; 1061
Keski-Rahkonen Pekka	375	Kim Young Chang	766
Kessler Nikolas	341	Kim Young Eun	264
Keum Soon Oh	1110	Kim Young Hwan	1051; 1061; 1146
Khaled Abir	992	Kim Young Soo	1084
Khalikova Maria	644	Kim Young-Hoon	560
Khamis Mona	217	Kim Young-Kwan	1195
Khamitova Gulzhan	1096; 1206	King Adam	719
Kharybin Oleg	157; 178	King-Ahmad Amanda	64
Khatib Soliman	1140	Kiontke Andreas	452
Khattar Rikkita	81	Kirchkeszner Csaba	982
Khomenko Luliia	34; 136; 370; 873	Kiriazis Alexandros	814
Khoo Amanda	436	Kiriş İrem	414; 525
Khoury Spiro	1128	Kirk Ansgar T.	943; 1157; 1176
Kiaei Mahmoud	384	Kirk Jayne	573; 861; 1217
Kiefer Karin	184	Kirkpatrick Joanna M	426
Kieffer Martin	733	Kishi Eri	1004
Kiehl Douglas E	823	Kishida Marina	510
Kiehl Tomotheus	1025	Kislinger Thomas	436
Kiesilä Anniina	1196	Kitata Reta Birhanu	436
Kihara Keishi	909	Kitayama Tetsuya	623
Kijewska Monika	165	Kitchen Simon	552
Kikuchi Yo	693	Kitov Pavel	168
Kilar Aniko	1121	Kitova Elena	168
Kilar Ferenc	1121	Kiviranta Hannu	333; 1063
Kilgour David P.A.	159; 631	Kiyonami Reiko	961; 1145
Kim Bo Kyung	714; 718	Kiyota Taira	490
Kim Bum Jin	659	Kjær Christina	309
Kim Byungjoo	1239	Kjeldsen Frank	90; 248; 280; 405; 884
Kim Chanmi	382	Klánová Jana	38; 432; 587; 616
Kim Dong-Hyun	496; 759; 760; 1152	Klassen John	168
Kim Eun Jeong	981	Kleberg Karen	906
Kim Eunjin	1197	Klee Sonja	580
Kim Geumsoog	1153	Klein Bruna	996
Kim Hugh I.	784	Klein Christian	604; 1182
Kim Hyo-Jin	262; 416; 717; 721	Klement Ales	1089
Kim Hyun Sik	332; 464; 466	Klieber Martin	368

Klijnstra Mirjam D.	1008	Kosgei Gilbert	303
Klodnitsky Helen	1138	Kosinsky Thomas	275; 649
Klont Frank	265	Kosior Piotr	1025
Klucken Christian	1082	Kosjek Tina	527
Kluijtmans Leo	725	Koske Daniel	143
Klute Felix David	811	Koster Sander	977
Kniemeyer Olaf	415	Kösters Steffen	266
Knížek Antonín	944	Kostiainen Risto	444; 633; 770
Knorr Arno	602	Kostmann Susann	795
Knotts Tom	604	Kostyukevich Yury	178; 468
Knox Ruth	172	Kosyakov Dmitry	1064
Knudsen Kristina B.	248	Kotiahho Tapio	444; 770; 814
Ko Jung Ah	1207	Kotimoole Chinmaya Narayana	603
Koal Therese	755; 1141	Koulman Albert	391
Kobayashi Hironori	381	Koundouros Nikolaos	193
Kobeissy Firas	148	Kovac Andrej	557
Koch Heiner	275; 427; 649; 1188	Kovacech Branislav	557
Koch Scarlet	275; 427; 649	Kovalev Sergey	591
Kochkin Dmitry	819	Kowalska Joanna	245
Kocsis Bela	1121	Kozanli Murat	592
Kocylowski Maciej	401	Kozhinov Anton N.	159; 631; 958; 962
Kodesova Radka	1089	Kožich Viktor	697
Koehler Karen	272	Kozlov Boris	593
Koehler Peter	410	Kragl Udo	1039
Kofoed-Sørensen Vivi	1062; 1193	Krah Philipp	635
Kohler Malcolm	216; 729	Krampera Mauro	288
Köhrle Josef	387	Kranjc Simona	527
Kokesch-Himmelreich Julia	910	Kratzke Manuel	755
Kolkman Annemieke	1050; 1094	Krause Michael	1188
Koller Dóra	526; 550	Krawczyk Bartłomiej	1172
Kollewe Astrid	401	Krawczyńska Anna	720
Köllner Bernd	267	Krechmer Jordan E.	580
Komarov Alexander	987; 1010	Kreutzmann Arne	224
Komatsuzaki Ryo	823	Krijt Jakub	697
Kompauer Mario	922	Krimm Pauline	876
Kondyli Aikaterini	589	Krisp Christoph	267
Konieczna Lucyna	380; 542; 720	Krizova Vera	110
Konijnenberg Albert	781; 799	Krmela Ales	983
Kononikhin Alexey	178; 289; 292; 459; 651	Kroeger Nils	92
Kooij Gijs	904	Kroesser Dennis	807; 867; 874
Kooij Pascal	1094	Kroniger Tobias	267
Kooijman Pieter C.	631	Kruathongsri Prapaijit	362
Kooistra Jeroen	439	Krüger Thomas	415
Kopaev Igor	590	Kruppa Gary	80; 427; 649
Koponen Jani	333; 1063	Kruth Claudia	457
Kopra Jaakko	633	Kruve Anneli	815; 832; 1218
Kopysov Vladimir N.	308; 309	Kruve-Viil Anneli	568; 1196
Korf Ansgar	1129	Kryštof Vladimír	544
Korfmacher Walter	130	Krzan Mojca	527
Kornilova Anna	467; 494	Kuballa Juergen	135
Korzeniowska Barbara	884	Kubeš Martin	728; 733
Kosek Eva	258; 283	Kubo Ayumi	1228

Kučera Lukáš	1240
Kucera Radek	385
Kucsma Nora	740
Kuepfer Lars	562
Kukačka Zdeněk	536; 541
Kuki Ákos	59; 470; 477; 503
Kuleff Alexander	109
Kulkarni Purva	626
Kultima Kim	1211
Kumano Shun	494
Kumar Ananthanarayanan	781
Kumar Anoop	994
Kumar Eva	1063
Kumar Seema	540
Kunath Tilo	745
Kuno Hideyuki	756
Kunz Laura	562
Kupcik Rudolf	425
Kupiec Monika	524
Kursinszki László	1038
Kurulugama Ruwan	1182
Kurumizaka Hitoshi	1033
Kurz Philipp	107
Kuuranne Tiia	140; 889
Kuwabara Hideya	1156
Kuzan-Fischer Claudia	191; 490
Kuzma Marek	1236
Kwan R.	429
Kwiatkowski M. D.	650
Kwiatkowski Marcel	334
Kwon Min-Hee	718
Kyung Cheol Yi	1084

L

La Rocca Cinzia	1053; 1117
Labouré Hélène	33
Laboureur Laurent	1144
Labuda Aleksander	497
Lackmann Jan-Wilm	636; 653
Lacko Michal	182; 944
Lacmanová Tereza	348
Lacoursiere Jean	889; 990
Laczko Endre	706
Ladocsi Lilla	1038
Ladumor Mayurbhai Kathadbhai	603
Lafiosca Maria Cristina	704
Lago Alexandre	307
Lai Bo	1060
Lai Cheuk-Kuen	129
Lalleman V.	556
Lam K. H. Brian	129
Lam Pui Yiu Yuko	75; 91; 158; 959; 1230
Lamb Aaron	974
Lamberti Cristina	1100
Lamkadam Younes	439
Lammert Stephen A.	114
Lammert Tobias	901
Lamonaca Fabio	831
Lamont Lieke	629
Lamoureux Catherine	880
Lancaster Cady	96
Landolfo Santo	16
Lanekoff Ingela	149
Lange Oliver	224
Lange Tobias	260
Langeland Knudsen Jeppe	309
Langley G. John	335
Langridge James	1190
Langsdorf Markus	755
Lantz Carter	155
Laouby Zahra	15
Laporte Delphine	1041
Lappo Nikita	1220
Lara Francisco J.	1165
Laramee Brittany	471; 495
Larcher Roberto . 161; 177; 235; 857; 1011; 1018	
Largo Antonio	112
Larraillet Vincent	645
Larsen Martin Røssel	436; 442
Larsson Anders	1211
Laserna Anna Karen Carrasco	394; 617
Laskin Julia	98
Latkin Tomas	1064
Lattanzi Claudia	666
Lau Adam	413
Lau Justin Kai-Chi	129
Laub Annegret	594
Laurent Sophie	1164
Lavanant Hélène	200
Lavold Thorleif	398
Law Richard	973
Lawal Remi	476
Lay Jaxkson Jr.	488
Lazarev Alexander	28; 536
Le Bizec Bruno	1165; 1166
Le Blanc Yves J. C.	806; 1173
Le Minh Uyen Thi	911
Le Quéré Jean-Luc	33
Le Thao Nhi	179; 1047
Le X. Chris	233
Lebedev Albert ... 311; 575; 591; 818; 1064; 1072	
LeBlanc Luc M.	684
Lech Katarzyna	1025
Lee Byoung Seob	332

Lee Chi Gyu	1065	Leogrande Patrizia	142
Lee Choong Sik	941	Leona Marco	7
Lee Choong-Min	560; 561	Lepine Franck	109; 298; 1180
Lee Dae Young	766	Leporati Marta	984; 1104
Lee Daeyoung	1153	Lerch Markus	1039
Lee Duk-Hee	718	Lesage Denis	99; 940
Lee Gwnagbin	1131	Leslie Elaine	233
Lee H. Y.	561	Lesniak Lukas	690
Lee Hangyeore	293	Leung Donald Y. M.	1115
Lee Holly	1000	Levent Serkan	821; 822; 825; 832; 833
Lee Hong Hee	1239	Levernæs Maren C. S.	378
Lee Hui-Ling	1067	Levi Mikaël	381
Lee Hyeyoung	1239	Lewandowska Aleksandra	417
Lee J.	561	Lewis Kathleen	413
Lee Jae Won	766; 1153	Lewis Michael	716
Lee Ji Hyun	597	Lewis Robert	574
Lee Jiyeong	416; 721	L'Hermite Jean-Marc	688
Lee John	604	Li Ailin	114
Lee Jong Cheol	1130; 1131	Li Enyou	32
Lee Juhyun H.	1066	Li Fong Yau Sam	617
Lee Kyoung-Seok	1109	Li Frederick	471; 495
Lee Maw-Rong	1113	Li Fu-An	646
Lee Sang-Won	293; 1235	Li Haiyang	32; 52
Lee Seok-Yong	560; 561	Li Huilin	155
Lee Seon eui	828	Li J.	429
Lee Seung-Ho	988	Li Jun	1060
Lee Sooyeon	372	Li Liang	748
Lee Sun Young	465; 642; 910	Li Meng	75; 158; 959
Lee Sung Bae	911	Li Sam F. Y.	394
Lee Tae Geol	911	Li Xiang	337
Lee Tung-Chun	208	Li Xianjiang	25
Lee Yong-Moon	1042	Li Xin	331; 844; 1132
Lee Yoo-Jin	262; 416; 717; 721	Li Xing-Fang	233; 1068
Lee Yoo-Rim	416	Li Yihan	535
Lee Youngseob	1153	Li Yiwen	337
Lee You-Rim	262; 717; 721	Li Zhiming	357; 359; 607; 1085; 1215
Lee Yun Jeong	561	Li Zhixiong	168
Leeming Michael	221	Lia Xiang	1060
Legros Véronique	499	Liang Di	389
Lehmann Fredrik	523	Liang Paul	471; 495
Lei Hua	32	Liang Tao	1123
Leito Ivo	351	Liang Yuxue	89
Lejeune-Dodge Aurore	1186	Liao Ruoying	337
Lemaire Joel	112; 208	Libralesso Elisa	831
Lemaur Vincent	1164	Lictevout Jean-Christophe	350
Lemaux Peggy	658	Lid Marthe	845
Lemos Fumiere Mayara	1103	Lieblein Tobias	790; 1169; 1187
Lemr Karel	491; 969; 989; 1168	Liew Chia Yen	1021
Lenčo Juraj	336	Liigand Jaanus	351; 815; 832; 1218
Lengyel Josef	301	Liigand Piia	351; 832; 1218
Leni Davide	271	Lim Dong-Kwon	465
Lento Cristina	789	Lim Hee Jin	465

Lim Heung Bin	382	Löfgren Lars	769
Lim Hui Fang	394	Loftus Neil	269; 615
Lim JongChoo	502	Logrieco Antonio F.	1006
Lim Moo-Song	981	Lok Natalie	287
Lim Myung Chul	1109	Lombard-Banek Camille	123
Lim Sang Ho	1065	Lomonaco Tommaso	37
Lim Youngran	1109	Long Michael C.	321
Lima Diogo Borges	775	Longo Edoardo	188
L'Imperio Vincenzo	920	Longo Valentina	366
Lin Cheng	683	Lonzarich Valentina	873
Lin Erika	726	Loo Joseph A.	155; 472
Lin Huan	1122	Loo Rachel Ogorzalek	472
Lin Linda	823	Looby Nikita	65; 383; 528
Lin Shanhua	788	Loos Martin	1049
Lin Yu	1231	Lopata Andreas	812
Lindblad Ulrik	634	Lopes Barbosa Ingrid	884
Linde Cecilia	769	Lopes Norberto Peporine	735
Lindner Herbert	209; 647	Lopez Felipe	515; 580
Ling Stephanie	147	López- González Ángeles	737
Linington Roger G.	1160	Lopez Nathan	1182
Linscheid Michael W.	942	López Ruíz Rosalía	1081
Linseele Veerle	1242	Lopez-Avila Viorica	619
Liou Saou-Hsing	1067	Lopez-Ferrer Daniel	137; 269
Lioznov Anton	157	Lorenz Christin	436
Liscio Camilla	745	Lorenz Matthias	497
Lisignoli Gina	253	Lorenz Sybille	434
Liskova Ruzena	10	Loriot Vincent	109
Lisowski Marek	1202	Losito Ilario	868
Liss Joanna	417	Lößner Christopher	312; 529; 647
List Katharina	618; 656	Louarn Essyllt	208; 1087
Little James L.	1046	Loula Martin	1070; 1198
Liu Baohong	624; 758; 1231	Loup Benoit	895
Liu Chang	677	Lourenço Célia	35; 367
Liu Huwei	25; 450; 965	Louzan Anna Letícia	189
Liu Jingyi	829	Lovrics Anna	740
Liu Mine-Yine	1133	Lowe Bethany	153
Liu Qingqing	233	Loyd Bill	806
Liu Shuang	791	Lozano Julian	1043
Liu Siqi	436	Lu Guohui	1071
Liu Suya	614	Lu Hui	793
Liu Wei	1069	Lu Xiufen	233
Liu Yu-Ching	637	Lubben Anneke	539
Livet Sandrine	440	Lubeck Markus	275; 427; 649; 1188
Liyanage Rohanna	488	Lucatello Lorena	975
Ljung Karin	728	Lukassen Marie V.	654
Llorca Marta	60	Łukaszuk Krzysztof	257; 417
Llorente-Onaindia Jone	554	Luna Baron Ana-Laura	208
Locachevic Gisele Aparecida	1137	Lunardi Andrea	738
Locatelli Marcelo	531	Lunkenbein Thomas	1082
Lock Stephen	312; 529; 647; 698	Lunn John E.	692
Lockhauserbäumer Julia	805	Luo Xian	748
Lockyer Nicholas	1134	Luo Xin	620; 621; 622; 1009

Lupu Loredana 28; 281; 536; 541; 807
 Luterova Andrea 1213
 Lutz Jean-François 58; 674; 1161
 Luu Betty 490
 Lynch Alice 964
 Lynch Iseult 209
 Lyskova Pavlina 379

M

Ma Qiang 55; 206
 Ma Shuguang 701
 Maagdenberg Arn 915
 Máčalová Eva 969
 Macaluo Veronica 128
 Maccelli Alessandro 939; 1105
 Mach Mateusz 830
 Mach Phillip 473; 581; 883
 Machálková Markéta 912
 MacIntyre Sam 1224
 Mackiewicz Elzbieta 664; 673
 Macklin Andrew 436
 Mackulak Tomas 605
 Macrì Antonella 1242
 Macur Katarzyna 257; 417
 Madden John E 602
 Madeira Jean-Paul 654
 Madera Mary 658
 Madsen Jeppe Buur 169
 Maeda Hiroto 277
 Maeser Stefan 536; 541
 Maestri Claudio 1034
 Magalhães Kelly 1139
 Maghari Shokoofeh 1199
 Magi Emanuele 313
 Magne Lionel 1087
 Magni Fulvio ... 253; 271; 273; 916; 920; 921; 923
 Magnússon Runa 97
 Magréault Sophie 880
 Magyar Anna 569
 Maharajh Rob 134
 Mahnke Lisa K. 106
 Maitre Philippe 306; 810
 Mak Tytus D. 167
 Makarov Alexander 156; 224; 590; 781
 Maker Garth 719
 Makino Elina 1138
 Malásková Michaela 368
 Maldaner Adriano Otávio 124
 Maldini Mariateresa 723; 724; 726
 Malhan Neha 658
 Malkawi Abeer 812
 Mallagaray Alvaro 805

Mallah Khalil 148
 Mallard Gary W. 1046
 Malutin Alexandr 324
 Malysheva Svetlana 985
 Managh Amy 929
 Manega Chiara Maria 753
 Maneta-Stavarakaki Stefania 474
 Manfredi Marcello 16; 273; 288; 327;
 418; 763; 1100
 Mangone Annarosa 1203
 Mangraviti Domenica 859
 Mann Graham 399; 400
 Mann Matthias 286; 767; 1188
 Mannerkorpi Kaisa 258; 283
 Mannina Luisa 1105
 Mannocci Francesco 963
 Mantovani Cynthia 189
 Manzano Sánchez Lorena 855
 Marabese Mirko 699
 Maranghi Francesca 1053; 1117
 Marangon Elena 513
 Marcelo Paulo 885
 Marchante-Gayón Juan Manuel 419; 422
 Marchetti-Deschmann Martina 913
 Marchis Daniela 984; 1104
 Marcillo Andrea 337
 Marciniak Alexandre 109
 Marcoux Judith 402
 Marcus Kathrin 266
 Marengo Emilio 16; 273; 288; 327;
 418; 763; 1100
 Maresch Daniel 1020; 1023
 Mariani Alessandro 1076
 Marín Sáez Jesús 836; 855; 979
 Marincas Olivian 1074
 Marini Federico 142
 Marino Gennaro 227
 Märk Lukas 478; 1080
 Markey Sanford 89; 1046; 1232
 Marlton Samuel 21
 Marquis Virginie 1006
 Marshall Alan G. 78; 201
 Marshall David L. 125; 192; 687
 Marshall John 147; 924
 Martens Jonathan 19; 294; 295; 725; 938
 Marti Matteo 881
 Martin Janosch 790; 1169; 1187
 Martin Jean-François 978
 Martinez Alexandre 99
 Martinez Jairo Rene 1043
 Martínez Vidal José Luis 1081
 Martínez-Jarquín Sandra 96; 489

Martinez-Lozano Sinues Pablo	216; 729	McArdell Christa	601
Martin-Garcia Beatriz	1237	McCall Holly	660
Martinotti Simona	273	McCort Gary	556
Martins Aline	1139	McCullagh James	791
Martins-Greens Manuela	751	McCullagh Mike	163; 1171; 1190; 1212
Marton Andras	875	McDonnell Liam	146; 285; 423; 679
Marton Marian	605	McDougall Danielle J	1134
Maruyama Hitoshi	494	McDowell Douglas	393
Marziali Sara	986	McDowell Melanie	116
Marzocchi Silvia	986	McGregor Laura	339; 875
Masaki Noritaka	475	McIntyre Will	187
Mascini Marcello	860	McKenna Amy M.	201
Masiello Lucia Anna	970	McKenzie James	898; 908
Masilamani Jeyakumar	617	McKinnon Benjamin	21
Masliuk Liudmyla	1082	McLeod Malcolm	813
Masnikosa Romana	839	McLuckey Scott A.	150
Massafra Stefania	984	McMahon Adam	900; 1134
Masselon Christophe	117	McManus Donald	389
Masson Elodie A. Y.	1128	McMartin Dena	1058
Masten Jasmina	420; 853	McMillan Daniel	595; 1073
Mastrorocco Fabrizio	1203	McPhaul Michael	285
Masuda Keiko	263	Mederos América	345
Masuro Domenico	705	Medina Mireia	847
Matabosch Xavier	554	Medvecka Erika	605
Mateus Julián A	554	Mehmood Shahid	116; 783
Mathiron David	834	Meid Andreas D.	518
Mätlik Kert	444	Meier Florian	286; 767; 1188
Matfoka Mikołaj	830	Meireles Maria Angela A.	1036
Matsui Takuya	756	Meitei Ningombam Sanjib	1142
Matsumoto Takuya	54; 330	Mejía Ospino Enrique	202
Matsuo Jiro	493	Melani Rafael	956
Matsuzaka Takashi	1149	Melchiorre Chiara	1032; 1244
Mattern Derek J.	337; 415; 1091	Melichar Bohuslav	1126
Mattivi Fulvio	705; 738	Mellors J. Scott	352; 960
Mattoli Luisa	571; 1037; 1099; 1242	Melnikova Nataliya	928
Matúšková Miroslava	545	Meltretter Jasmin	979
Maunit Benoît	179; 1047	Melymuk Lisa	587; 616
Mauracher Andreas	944	Memboeuf Antony	686
Maxouti Katerina L.	945	Memo Maurizio	784
May Caroline	266	Mendez Andrea	208
Mayer Paul	153; 690	Menegon Jair	302
Mayhew Chris	887	Menezes Marcel	826
Mayhew Christopher	368	Menezes Riya Christina	434
Mayo Barrallo Juan Carlos	737	Mengerink Ynze	57
Mayr Christine M.	1106	Menghi Leonardo	370
Mayrhofer Corina	441	Menicatti Marta	840; 841; 846
Mazur Dmitrii	818; 1072	Meniche Xavier	716
Mazza Mariafrancesca	253; 920; 921	Menin Laure	159; 1226
Mazzarino Monica	877; 881	Mercati Valentino	1242
Mba Cecile	550	Meriaux Celine	15
Mc Cann Andréa	338	Merkel Dietrich	413; 726; 1170
McAlister Graeme	428; 731; 823; 961; 1184	Merkyte Vakare	188

Merlanti Roberta	975	Mocsai Reka Tunde	1024
Messinger Jeffrey	917	Modenbach Jana Marielle	1039
Mestdagh Hélène	1087	Modzel Maciej	1202
Mestek Oto	1070; 1198	Moeller Alexander	454
Metalnikov Pavel	987; 1010	Moellers Rudolf	628
Meurs Joris	496	Mohren Ronny	420
Meyer Bernd	1031	Mohtashemi Iman	731; 823
Meyer Sven	767; 1142; 1177	Moilanen Jani	1196
Michael Chalmers	171	Mok Hyejung	1197
Michel Alain	864	Mokh Samia	531
Michopoulos Filippos	769	Moldovan Zaharie	1073; 1074; 1241
Micic Srdjan Johnny	454	Molinaro Francesco	540
Mielcarek Andreas	604	Molinié Roland	834
Migas Lukasz	173; 1219	Mollahosseini Afsaneh	1194
Miggiano Riccardo	514	Møller-Jensen Jakob	442
Miglio Roberta	1147	Molloy Billy	278
Mikhailov Victor	791	Molloy Mark P.	267
Mikhalychev Alexander	1220	Molnár Adrienn	340
Mikkilä Jyri	461	Monaci Linda	868
Mikoviny Tomas	1080	Monastyrskiy Mikhail	590
Milani Hosseini Mohammad-Reza	1194	Mondello Luigi	237; 859
Milella Luigi	238	Monfort Jordi	554
Milentyev Alexander	92	Monfort Nuria	554
Milet Anne	940	Monk Raymond	387
Mileto Simona	230	Montagna Daniela Domenica	366
Miller Ken	435	Monteau Fabrice	1165; 1166
Milne Joy	745	Monteiro Olinda	1054; 1075
Min Ayoung	597	Montemurro Nicola	1078
Min Hyung Sik	1109	Monti Maria Chiara	408; 532
Minameyer Martin B.	945; 1200	Moody Sally	71; 123
Minassi Alberto	532	Moon Dae Won	465; 910
Miolo Gianmaria	364; 703	Moon Jeong Hee	911
Mirabelli Mario Francesco	475; 768	Moon MyeongHee	1130; 1131
Miralles-Marco Ana	587	Mooney Mark H.	144
Miranda Renata	884	Moore Jessica	122
Miraval Tommaso	530	Moorthy Arun	95; 1221
Mirokhin Yuri	89; 1232	Moraes Luiz Alberto Beraldo de	1137
MIronova Polina	250	Morag Ely	798
Mishima Yuji	353	Morel Nathalie	256
Mitchell Todd W.	687	Moreno-Pedraza Abigail	489; 630
Mitro Nico	219; 694; 734	Moretti Simone	995; 1056
Mitrofanov Elena	680	Moretto-Capelle Patrick	688
Mitsui Shinichi	710	Morgan Tomos	158; 964
Mittal Parul	121	Morgner Nina	77; 241; 462; 685; 786; 790; 799; 1169; 1187
Miura Daisuke	930	Morineau Loïs	906
Miyahara Miho	276	Morishima Ken	802
Miyazaki Masaru	276	Moritz Franco	80; 596
Mizuno Hiroki	494	Morozova Ksenia	362
Mladěnka Přemysl	702	Morrice Nick	45; 413
Mladenovic Katarina	1011	Mortishire-Smith Russell	1171; 1172; 1190; 1212; 1217
Mochalova Elena	882; 892		
Mochalski Pawel	368		

Mosely Jackie 358
 Mosen Peter 421
 Moser Sergio 857; 1011
 Moshkovskii Sergei 280
 Mosler Dominik 491
 Moszczyński-Pętkowski R. 830
 Motamedchaboki Khatereh 726; 741
 Motoyama Akira 909
 Mouchet Elizabeth 492
 Moulard Yves 893; 895
 Mourad Daniel 224
 Moyano Encarnación 481; 576; 847
 Mozzicafreddo Matteo 532
 Mroczek Tomasz 414; 525
 Muccilli Vera 866
 Mueller Simona 454
 Mulder Inge 915
 Mulinacci Nadia 705
 Mullen Christopher 956
 Müller Adrian 184
 Müller Markus 1080
 Müller Max 922
 Müller Oliver J. 518
 Müller Sebastian 250
 Müller Tobias 415
 Müller Volker 799
 Mullin Lauren 719
 Mun Dong-Gi 1235
 Mun Sora 262; 416; 717; 721
 Munem Marwa 934
 Munier S. 864
 Munoz David 490
 Muñoz de la Peña Arsenio 341; 854; 1237
 Munshi Musleh 295
 Muntel Jan 45; 250; 426
 Murakami Taro 1004
 Murphy Keeley 1138
 Murphy M. 820
 Murray David 492
 Murray Kermit K. 86; 476
 Murta Teresa 147; 924; 1222
 Muruve Daniel 921
 Muselli Alain 361
 Muskat Tassilo 680
 Musselman Brian 471; 495
 Mutschlechner Paul 1080
 Myers Gabriel 658
 Myler Peter 956
 Myśliwiec Małgorzata 380
 Myung Seung-Woon 988

N

Na Seungjin 1223
 Nagadoi Aritaka 1033
 Nagano Hisashi 494
 Nagano Maiko 510
 Nagatomo Kenji 1228
 Nagayama Satoshi 263
 Nagornov Konstantin O. 159; 631; 958; 962
 Nagy Lajos 59; 470; 477; 503
 Nagy Tibor 59; 470; 477; 503; 570
 Najmanová Iveta 702
 Najumudeen Arafath K. 147
 Nakamura Sayaka 504
 Nakamura Takemichi 342
 Nakano Yosuke 727
 Nakata Keisuke 564
 Nakayama Hiroshi 9
 Nakayama Kenji 1132
 Naldi Marina 281; 343; 562
 Ñañez Nathalie 1043
 Napoletano Sabino 669; 670
 Nappi Raffaella 972
 Nardin Tiziana 161; 177; 1018
 Narmandakh Mijid 59
 Narumi Ryohei 263
 Nattkemper Tim W. 901
 Naudé Yvette 29; 344
 Navarini Luciano 865; 873; 1096
 Navarrete Perea Jose 428
 Navarro Pedro 435; 436
 Navarro Suarez Laura 711
 Navrátilová Jarmila 912
 Nawrocki Arkadiusz 442
 Nebbia Sterfano 1100
 Nebl Thomas 1224
 Neittaanmäki Noora 934
 Nekouei Mina 384
 Nemat Reza 187
 Nemes Peter 71; 123
 Neri Bruno 613
 Neri Dario 520
 Neta Pedatsur 757
 Neto Francisco 884
 Netzer Ravit 780
 Neugebauer Thomas 681
 Neumann Anika 82
 Neumann Steffen 709; 713; 1211
 Neuweger Heiko 341
 Neveu John 69

New Leesun	1107	Novara Maria Eugenia	1204
Newman Carla	70	Nowak Nora	729
Newton Kenneth	1182	Nshanian Michael	155; 472
Ng Bao Hui	394; 617	Ntai Ioanna	730; 731; 746
Nguyen An	190	Ntasi Georgia	227; 1244
Nguyen Binh Thanh	534	Núñez Nerea	856
Ni Chi-Kung	1021	Núñez Oscar	139; 585; 850; 856; 1079
Ni Guoying	389	Nussbaumer-Ochsner Yvonne	216
Niazi Sadia	963	Nuti Francesca	165
Nicklay Joshua J.	435	Nyiri Zoltán	982
Nicol Edith	631	O	
Nicolás Carcelén Jesús	419; 422	O'Connor Peter B.	75; 91; 158; 959; 964; 1230
Nicoli Raul	889	O'Sullivan Gwen	1058
Nie Honggang	25	O'Brien-Johnson Reid	472
Niehaus Karsten	901	O'Doherty John V.	186; 1097
Niehoff Ann-Christin	914	Oertel Reinhard	242
Niehuis Ewald	628	Ogawa Osamu	1132
Nielen Michel W.F.	127; 160	Ogrinc Nives	234; 420; 848; 853; 862
Niemikoski Hanna	143	Ogrinc Potocnik Nina	915
Nieto Sofia	619	Oh Myung Jin	535
Nikolaev Evgeny	40; 157; 178; 289; 292; 468; 651; 928; 1151; 1209; 1227	Ohk Boram	714; 718
Nikula Chelsea J.	147; 924; 1222	Öhrfelt Annika	249
Nikulin Vyacheslav	987	Ohtsuka Masayuki	276
Niles Sydney F.	78; 201	Ojakivi Mari	815; 1218
Nilsson Ralf	769	Okada Yutaka	1201
Nilsson Sofia M. E.	814	Okazaki Shigetoshi	475
Nishimoto Haruka	369	Okimi Takuya	693
Nishimoto Koshiro	919	Okońska Magdalena	380
Nishimura Mai	710	Okuda Koji	321; 1126
Nishimura Yoshifumi	1033	Oldekop Maarja-Liisa	351
Nishiwaki Yoshinori	710	Ołdziej Stanisław	257; 417
Niu Zengyuan	620; 621; 622; 1009	Olivaro Cristina	345; 346
No KiMi	981	Oliveira Luan Felipe C.	689
Nöbauer Katharina	441	Oliveira Maria da Conceição	1075
Nocentini Alessio	841	Olivenza-Leon David	368; 887
Nogueira Gislaine C.	1036	Olkowicz Seweryn	38
Noij Theo	439	Olmo-García Lucía	341; 854
Nojiri Tatsuo	494	Olsen Jesper Velgaard	88
Noll Thomas	404	Olsen John Elmerdahl	442
Nomura Fumio	276	Olson Loren	754
Nonami Hiroshi	564	Olszewski Maureen	1138
Norais Nathalie	1034	Oltmanns Sven	204
Nørgaard Asger W.	1062; 1193	Omer H.	654
Nortcliffe Chris	535; 543; 648; 785; 1173	Oncak Milan	301
Nothias Louis Felix	582; 670	Ončák Milan	20; 104; 107; 304
Novaes Fabio J. M.	1114	Ondruchová Jana	1168
Novak Jiri	1236	Onjiko Rosemary	123
Novák Ondřej	728; 733; 932	Onzo Alberto	704; 732
Novak Petr	10; 385	Oomens Jos	19; 129; 294; 295; 725; 938
Nováková Lucie	348; 702	Opalka Lukas	38
Novakovskiy Roman	928	Operti Lorenza	863
		Oppenheimer Jim	619

Oppermann Christina	1039	Palfner Götz	594
Oppermann Madalina	792	Pallecchi Marco	395
Organtini Kari	861	Palm Martin	933
Ori Alessandro	426	Palmer Martin	1190
Ori Carlo	39	Palmieri Marta	763
Ornoch Anna	1025	Palmisano Francesco	316; 578
Orsini Paolo	548	Palombo Paolo	1003
Orsomando Giuseppe	512	Palsson Runolfur	396
Ortiz Alexia	644	Panina Yulia	924
Ortiz Daniel	1226	Panitchpakdi Morgan	582
Ortori Catherine	743	Pannequin Anaïs	361
Osakabe Akihisa	1033	Paola Signorelli	374
Osawa Rodrigo	1054; 1075	Paoella Francesca	253
Oschwald Johannes	945	Paoli John	934
Oshikawa Tomonori	930	Paolini Mauro	857
Osmaniye Derya	833	Papan Cyrus	723; 724; 726; 1120; 1170
Oster Andréia Hansen	1111	Papanastasiou Dimitrios	682
Otake Yukiko	369; 1192	Papini Anna Maria	165; 536
Oteri Marianna	237	Papp Dávid	650
Othman Alaa	764	Pappin Darryl J.	446
Otsuka Yoichi	54; 330	Parczewski Andrzej	664; 673
Otto Andreas	267	Parent Annick	556
Otzen Daniel	169	Pařízková Barbora	733
Oulie Inger	386	Park Arum	262; 416; 717; 721
Ovando-Vázquez Cesaré	489; 630	Park Gun Wook	1217
Ovchinnikova Olga	497	Park Hyoung-Joon	597
Oyler Benjamin	1123; 1135	Park Jeong Ju	941
Ozaki Asako	1004	Park Jin Young	659
Ozalp Burhan	347	Park Jinkyu	1065
Ozan S. D.	1045	Park Ji-Won	465
Ozcan Saniye	832	Park Jong-Ho	1016
Ozdinler Pembe Hande	291	Park Kihong	1061
Özkay Ümide Demir	821; 825	Park Moon-Hee	1051; 1061
Özkay Yusuf	822; 825; 833	Park Seongsoo	597
Ozohanics Oliver	1028	Park Sung-Kwan	597
Ozturk Yusuf	821	Park Yong Chjun	1014
P			
Pačes Ondřej	354	Park Young Ja Hwang	1146
Pacholarz Kamila	67; 793	Parker Benjamin L.	436
Paciotti Roberto	810	Parker Robert	267
Paczkowski Dariusz	712	Paron Igor	286
Padalino Massimo	744	Paroni Rita Clara	374; 392; 739; 1138
Paek Eunok	1223	Paša-Tolić Ljiljana	81; 658
Paetz Christian	434	Pasca Marilina	418
Pagiotti Rita	1242	Pascal Quentin	458
Pagni Fabio	271; 916; 920	Pascale Raffaella	704; 732
Paine Martin	79	Pascher Tobias	301
Pala Anna Paola	1053	Pascovici Dana	267
Palacio Lozano Diana Catalina	202	Pasquali Marzia	36
Palarasah Yaseelan	90	Pasquatti Taynara	189
Palaric Cécile	834	Pasquier Stéphane	1087
		Passig Johannes	181; 457
		Passmore Lory A	781

Passoni Alice	1076	Pereira Henrique M. G.	1114
Pastorelli Roberta	699	Pereverzev Aleksandr Y.	308; 309
Pásztor Pavel	182	Perez Consuelo	1103
Patel Bhavin	410	Pérez Cova Miriam Carolina	1136
Patel Ekta	665	Pérez Sandra	1078
Patiny Luc	1226	Pérez-Mañá Clara	554
Paton Nicholas	394	Pérez-Mellor Ariel	48
Patriarca Claudia	1077	Perini Matteo	1017; 1018
Patsilnakos Alexandros	946	Perlman David	69
Patterson N. Heath	917	Permentier Hjalmar	794
Patterson Nathan	122	Permi Perttu	765
Pauers Martin	618; 656	Perminova Irina V	178
Pauk Volodymyr	989	Perret Alain	750
Paul P.	556	Personne Stephane	885
Pauling Josch	79	Persson Bengt	769
Paulo João	428	Peru Kerry	1058
Paulová Markéta	697	Perucchini Sara	1147
Paulus Aran	730	Peškova Karolína	697
Pavlatovská Barbora	912	Peter Ronja	663
Pavlík Jakub	348	Peterman Scott	435
Pawlaczuk Aleksandra	668; 888	Peters Ron	57
Pawlak Katarzyna	1025	Peters Ruud	211
Pawliszyn Janusz	23; 53; 65; 160; 383; 485; 486; 528; 992	Peters Thomas	805
Payen Didier	256	Peti Ana Paula Ferranti	1137
Peake David A.	1145	Petit Leseurre Lucie	179; 1047
Pearson Mackenzie J	1120	Petracca Roberto	1034
Pecci Alessandra	230	Petrik Milos	379
Pedemonte Nicoletta	403	Petrovics Noémi	982
Pedersen-Bjergaard Stig	845	Petry Stefan Matthias	549
Pedretti Silvia	694; 734	Pettit Michael	86
Pedrini Paolo	1015	Pezzati Paola	255
Pedro Liliana	565	Pfamatter Sibylle	44
Pedrotti Michele	34; 370	Pham Hai Tuan	755
Peez Oliver	77; 786	Pham-Tuan Hai	1141
Pekov Stanislav	40; 1151; 1227	Phan Toan Thang	617
Peleg Yoav	798	Phaneuf Clifford	1138
Pellarin Michel	688	Philip Robin	147; 924; 1222
Pellegrini Davide	423	Phillips Nancy	42
Pellegrini Sheila	881	Pianezze Silvia	1017; 1018
Pelupessy Philippe	966	Piao Haitao	1071
Peña-Herrera Juan Manuel	1078	Piatti Luca	1091
Pencik Ales	728	Picard Pierre	889; 990
Peng Hanyong	233	Picardo Massimo	585; 1079
Pengelley Stuart	559; 1158	Piccirillo Susanna	810
Penna Tatiana C.	302	Piccolo Stefano	12
Pennazio Enrico	984	Picenoni Renzo	768
Penoni Andrea	739	Pickford Russell	735
Pepi Federico	946; 947	Picotti Paola	17
Peporine Lopes Norberto	1160	Picron Jean-Francois	985
Perdones-Montero Alvaro	193; 474	Pieczkolan Jerzy	830
Perego Rita	548	Piel Felix	368; 478; 944; 1080
		Pierattini Erika Carla	1088

Piersanti Arianna	991; 1001; 1002; 1003;1108; 1234	Popov Mark	1064
Piersma Sander	436	Porchia Marina	105
Piestansky Juraj	557	Porro Benedetta	753
Piga Isabella	253; 271; 916	Porta Tiffany	460; 629
Pignatello Rosario	374	Portelius Erik	954
Pijnappel Matthijs	523	Porter Jacob	199
Pilard Serge	834; 885	Portero Erika	71
Pilařová Veronika	348; 702	Posocco Bianca	513
Pilkington Colin	844	Possi Alessandro	670
Pillai Manoj	994	Potapov Alexander	40; 1209; 1227
Pilon Alan Cesar	735	Poteshin Sergey	1174; 1179
Pimentel-Moral Sandra	1237	Potthoff Alexander	498
Pinkasová Renata	697	Poucheret Patrick	864
Pinto Gabriella	1032	Poulogiannis George	147; 193; 924
Pinto Provenzano Sara	366	Poulsen Nicole	92
Piper Thomas	349	Pouwels Simon D.	265
Pirkl Alexander	628	Poyer Salomé	58; 674; 1161
Pitonzo Rosa	1204	Prache Nolween	132
Pittalà Maria G. G.	657	Prada Fausto	1043
Pittenauer Ernst	100; 691	Prado Morgana Kelly Borges	1137
Pittia Paola	860	Prasad Satendra	428; 1184
Pitto Marina	273	Preisler Jan	505; 912
Piubelli Luciano	541	Prest Harry	619
Platania Chiara	374	Preud'homme Hugues	501
Pleiner Sina	272	Pries Antonia	404
Plet Benoit	580	Prieto Mari	137
Płóciennik Halina	1201	Pringle Steven	740; 875; 1219
Ploskonka Ann	473	Priola Emanuele	863
Pluhacek Tomas	379	Proczko-Stepaniak Maria	720
Plum Sarah	266	Progent Frederic	350
Plumb Robert	719	Proksch Roger	497
Plyushchenko Ivan	736	Prop Arnoud	915
Poad Berwyck	192	Pruška Adam	912
Poad Berwyck L. J.	125; 687	Przybylski Michael	28; 281; 536; 541; 807
Pogan Ronja	800	Pucci Piero	227
Pohl Kerstin	521; 535; 543; 648; 785; 1173	Pugh Jonathon	574
Pohnert Georg	907	Puignou Lluís	139; 850
Pöhö Päivi	770	Pukala Tara	121
Polasek Miroslav	110	Pulkrabova Jana	983
Polfer Nicolas	102	Punzi Angela	578
Poli S.	857	Purdy Georgiana	1150
Pollastro Federica	532; 1048	Purgel Mihály	59
Pollegioni Loredano	541	Putman Jonathan C.	78
Polyakova Olga	1072	Putz Marlen	349
Poma Genin Elisabetta	972	Puustinen Anne	669
Pompach Petr	385	Pyszka Magdalena	380
Ponzetto Federico	889		
Pop Nastasia	1241	Q	
Popolo Laura	281	Qiao Liang	624; 758; 1231
Popot Marie-Agnès	893; 894; 895	Qin Feng	337; 820; 967; 992
Popov Igor	40; 289; 292; 1151; 1209; 1227	Qin Qin	758
		Qing Feng	1091

Qu Lianyi 622
 Quack Thomas 922
 Quaglia Milena 445
 Quanico Jusal 15; 148
 Queguiner Laurence 538
 Querci Lorenzo 232
 Quercioli Massimo 255
 Quinn Ronald 956
 Quinton Loïc 338
 Quye Anita 228

R

Raab Shannon 195; 1035; 1175
 Rábai József 340; 570
 Rabezana Roberto 863
 Rack Jeffrey 303
 Radaszkiewicz Katarzyna Anna 912
 Radchenko Tatiana 1217
 Raddatz Christian-Robert 1176
 Radford Sheena 779
 Rados Edita 691
 Raether Oliver 649; 1188
 Raffaelli Andrea 1053; 1088
 Rafols Pere 751; 905
 Raftery Daniel 735
 Raggi Andrea 212; 213
 Raghuraman Bharath Kumar 424
 Ragno Rino 946
 Raguvaran Vanaja 979; 1000
 Rahmani Fereidoon 1194
 Raibaud Sophie 556
 Räikkönen Heikki 633
 Raimondo Francesca 273
 Rainville Paul 719
 Raja Erum 792
 Rajcic Boris 839
 Rak Alexey 556
 Rakov Dmitry 882; 892
 Ramanathan Dil 1011
 Ramanathan Ragu 64
 Ramirez Noelia 751
 Randak Tomas 453
 Randall Kristen 1138
 Rantakokko Panu 1063
 Ranzato Elia 273
 Rapinoja Marja-Leena 669
 Räsänen Riikka-Marjaana 633
 Rathmann Daniel 387
 Raupers Björn 479
 Raus Peter P.M. 269
 Rauschenbach Stephan 635
 Ravensborg Christian 286

Ray Andrew 358
 Ray Dilip K. 238
 Re Federica 253
 Re Nazzeno 810
 Rebane Riin 351
 Redman Erin 352; 960
 Reeksting Shaun 539
 Rehulka Pavel 425
 Rehulkova Helena 425
 Reich Fraser 887
 Reichel Christian 284
 Reid Gavin 43; 687
 Reiding Karli 166
 Reiner Eric 616
 Reinhard Jahn 795
 Reinhardt-Szyba Maria 156; 781
 Reiss Julius 635
 Reiter Lukas 45; 250; 426; 431
 Rektorisova Michaela 598
 Remita Hynd 208
 Remoroza Connie A. 167; 858
 Ren Hong 1167
 Renaud Tristan 834
 Renois-Predelus Gina 47; 1027
 Renslow Ryan 196
 Rentel Claus 243
 Repoux Marie 33
 Reubsæet Leon 378; 386
 Révész Ágnes 686; 1028
 Reyes-Garcés Nathaly 383
 Reynolds James C. 356; 773
 Rezaei Human 1163
 Rezende Claudia M. 771; 1114
 Ribeiro F. W. M. 300
 Riboni Nicolò 451
 Riccadonna Samantha 705
 Riccardi Sirtori Federico 540
 Ricci Andreina 946; 947
 Ricci Francesca 520
 Riccio Raffaele 408; 532
 Ricciutelli Massimo 24
 Richards Keith 387
 Richardson Susan D. 180
 Richer Stéphanie 1041
 Ricketts Christopher 1018
 Riedel Jens 484; 942
 Riera Antoni 1217
 Rigano Francesca 237; 859
 Rigaud Sébastien 834
 Righetti Laura 198; 1007
 Rilak Ana 839
 Rinaldi Francesca 541; 807

Ring Trial Consortium	80	Rossetti Fabrizio	188
Ringvold Amund	386	Rossi Franca	514
Risiglione Pierpaolo	657	Rossi Rosanna	995; 1056
Rissanen Matti P.	461	Rosta Edina	308
Riva Simona	540	Rosting Cecilie	378
Rivera Jeff Rivera	1000	Rosu Frédéric	200
Rizos Timolaos	518	Rosu Frédéric	245
Rizzo Jessica	374; 392; 739; 1138	Roszkowska Anna	65; 528; 542
Robert Stéphanie	733	Roth Adrian	562
Robert Vincent	99	Rousseau Alice	208
Roberts Gareth	875	Roussis Stilianos G.	243
Roberts Simon	1073	Rowland Steven M.	78; 201
Robinson Carol V.	73; 116; 781; 783	Rubert Josep	738
Robinson Kenneth N.	147; 1222	Rubino Federico Maria	374; 392; 739; 1138
Robitaille Aaron	433	Rudaz Serge	889
Robotti Elisa	16; 273; 288; 327; 418	Rudewicz Patrick	565
Rocchi Rachele	860	Rudolph Heather	742
Rodas Sánchez Laura Covadonga	737	Ruff Matthias	1049
Rodet Franck	15	Rüger Christopher	82
Rodgers Mary T.	8	Ruland Holger	1082
Rodgers Ryan P.	78; 201	Runolfsdottir Hrafnhildur L.	396
Rodin Igor	736; 819; 835	Ruoppolo Margherita	406
Rodrigues David A.	64	Ruppert Thomas	421
Rodrigues de Morais Damila	317	Rusche Hendrik	541
Rodrigues-Oliveira André F.	300; 302	Ruskic David	768
Rodríguez González Pablo	419	Russ Bill	619
Rodríguez Pilar	847	Russo Daniela	238
Rodríguez-González Pablo	422; 737	Russo Domenico	253
Roelen Bernard	563	Ruthenberg T.	820
Rogel Estrella	204	Rutigliano Barbara	972
Rogers John C.	433; 438	Rutka James	191; 490
Rogowski Jacek	664; 673	Rutkauskas Marius	801
Röhling Ulrich	87	Rutsch Moreno	601
Rohman Mattias	492	Ruzicka Viktor	385
Rohn Sascha	135; 641; 711	Ryaboshapkina Maria	769
Rohrer Jeffrey S	602	Ryan Daniel	122
Rohwer Egmont	29; 600	Rybinska Urszula	1012
Rokob Tibor Andras	686	Rychlik Michael	980
Rolando Christian	966	Ryumin Pavel	806
Rolli Enrico	1007		
Romanzin Claire	950	S	
Romão Wanderson	996	Sá Najla M. D.	771
Romero González Roberto	836; 855; 1081	Saaidi Pierre Loïc	750
Römpp Andreas	910	Sabatini Francesca	228
Rontogianni S.	408	Sabbah Hassan	678
Root Katharina	74	Sadler Peter J.	91; 1230
Root Yuriko	1011	Saffell John	35
Rosas-Román Ignacio	489; 630	Sage Ashley	994; 1107
Rosenfeld Mark	651	Sağlik Begüm Nurpelin	822; 833
Rosi Marzio	947; 951	Sagratini Gianni	24; 1096; 1102; 1206
Rosinha Gracia	189	Sahasrabuddhe Aniruddha	778
Rossell David	222	Sahin Cagla	169

Saikusa Kazumi	1033	Sarvin Boris	819; 835
Sainz Menéndez Rosa M.	737	Sarycheva Anastasia	1179
Saito Naoaki	809	Sato Futoshi	342
Sakai Mariko	623	Sato Hiroaki	500; 504; 506
Sakai Yoshiharu	263	Sato Hiroko	936
Sakakiyama Joji	328	Sato Kei	1167
Sakakura Motoshi	891; 1216	Sato Nobuhiro	802
Sakane Iwao	961	Satoh Takaya	500; 1228
Sakson Roman	421	Saudemont Philippe	458
Salanoubat Marcel	750	Saugy Martial	889
Salazar Gómez Jorge Iván	1082	Saurina Javier	139; 850; 856
Salehi Mehraveh	187	Savini Sara	212; 213
Saletti Rosaria	657; 866; 1112	Savio Valentina	970
Salomone Fabrizio	520	Sayers Rebekah	413
Salpin Jean-Yves	244	Scagliola Alessandra	273
Salsa Matteo	548	Scalbert Augustin	215
Salustriano Nathacha	996	Scampicchio Matteo	362
Saluti Giorgio	995; 1056	Scaramuzza Rosa Teresa	37
Salvatori Sandro	1242	Scarpone Rossana	669; 670
Salvitti Chiara	946; 947; 948; 951	Scarselli Maria	1034
Salzet Michel	15; 148; 458	Scavenius Carsten	654
Samarah Laith	81	Schäfer Mathias	787
Samguina Tatiana	591	Schäfer Ricarda A.	1200
Samino Sara	751	Schäffer Richárd	740; 875
Sammon Chris	665	Schalley Christoph A.	106; 568; 1196
Samonig Martin	730	Schar Manuel Y	615
Sanchez Ferran	521; 535; 543; 648; 785; 1173	Schatzmayer D.	968
Sánchez-Arévalo Carmen María	341; 854	Scheel Dierk	707; 709; 713
Sanda Miloslav	1026	Scheeren Simon Gereon	1162
Sander Julia	1101	Schefer Simone	641
Sander Peter	1177	Scherer Rodrigo	996
Sandor Viktor	1121	Scherf Katharina	410
Sandoval Wendy	776; 788	Scheuermann Jörg	520
Sandra Koen	644	Schey Kevin	917
Sansom Owen J.	147; 924	Schiller Sven Arne	1080
Santhanam Balaji	781	Schindler Baptiste	1027; 1180
Santorelli Lucia	273	Schindler Baptiste	47
Santoro Massimo	875	Schirinzi Gabriella Francesca	60
Santos Fabio	1139	Schlapbach Ralph	562
Santos Francisco Javier	481; 847	Schlemmer Stephan	299
Santos Vicente Francisco Javier	576	Schleuder Detlev	741; 997
Santovito Elisa	1006	Schlich Pascal	33
Saraji Mohammad	203	Schlögl Robert	1082
Sarbu Mirela	1030; 1178	Schlosser Gitta	340; 569; 570; 650
Sardela Vinicius F.	1114	Schlottmann Florian	1157
Sarg Bettina	647	Schlüter Hartmut	260; 334; 411; 612; 807; 867; 874
Sargenti Marta	1056	Schmid Robin	897; 1162
Saris Wim H. M.	250	Schmidt Carla	795; 801
Sarracino Dave	410	Schmidt Henning	113
Sarrazfadeh Mehrnaz	480; 494	Schmidt Jürgen	594
Sartorato Ambra	975	Schmidt Michaela	663
Sartore Luciana	253		

Schmit Pierre-Olivier	275	Sentellas Sonia	856
Schmitt-Kopplin Philippe	80; 596; 704; 732	Seo Jungju	1127
Schmitz Katja	28; 536	Seo Nari	535; 1027
Schneeberger Damaris	729	Seok Ae Eun	262; 416; 717; 721
Schneider Bernd	434	Seong Sook Jin	714; 718
Schnell Anne	816	Seong Yunseo	261
Schnitzler Jörg-Peter	596	Serafini Mauro	1098
Schoenmakers Peter	57	Serani Andrea	851
Schofield Christopher	791	Sergeeva Victoria	289
Schollee Jennifer	601	Sergi Manuel	582; 669; 670; 671; 860; 1098
Scholz Birte	452	Seró Raquel	481
Scholz Michael	303	Sersa Gregor	527
Schoorl Jorien	97	Shahhoseini Fereshteh	1057; 1083
Schorr Pascal	1181	Shakhmatov Dmitry	736
Schoutsen Frans	602	Shamshiddinova Maftuna	1042
Schrader Wolfgang	205; 589; 1191; 1229	Shan Paul	427
Schulte Uwe	401	Sharma Seema	701; 731; 823; 961
Schulz-Bull Detlef	181	Sharon Michal	115; 777; 780; 798
Schulzova Vera	983	Shaw Jared B.	81
Schüth Ferdi	205	Sheardova Katerina	38
Schwab Anatol	74	Shekargoftar Masoud	505
Schwaninger Markus	764	Shen Xiaopan	357; 359; 607
Schwartz-Zimmermann Heidi-Elisabeth	968	Shestivska Violetta	182
Schwarz Helmut	214	Shevchenko Andrej	92; 424
Schweiger Hufnagel Ulrike	748; 767; 1142	Shiea Jentaie	318; 918; 998
Schwepe Devin	428	Shieh Ying-Tzu	1133
Sciubba Fabio	669	Shiel Jonelle	311
Scolyer Richard A.	267	Shih Wen-Ting	1052
Scoppetta Fausto	1056	Shimabukuro Yuji	949
Scortichini Giampiero	1002; 1003	Shimano Hitoshi	1149
Scott-Boyer Marie-Pier	402	Shimbo Kazutaka	693
Scotti Gianmario	814	Shimizu Koji	1132
Scrucca Luca	1234	Shimizu Yoshihiro	263
Scrutton Nigel	1186	Shimma Shuichi	919
Scuderi Debora	112; 128; 244	Shin Chan Young	911
Scurr David	628	Shin Dong Woo	597
Scurti Elena	770	Shin Hyunsuk	275
Sears Vaughn	665	Shion Henry	573; 574
Sebastiani Luca	1088	Shiota Teruhisa	891; 1216
Sebela Marek	425	Shirshin Evgeny	178
Segura-Carretero Antonio	1237	Shishido Tânia Keiko	765
Seibold Max A	1115	Shizuma Motoshiro	1004
Seidel Ralf	801	Shokoufi Nader	365
Sekera Emily	742	Shon Hyun Kyong	911
Sekimoto Masayo	342	Shopova Iordana A.	415
Sekowski Jennifer	883	Shpigun Oleg	736
Selevsek Nathalie	562	Shurkhay Vsevolod	40; 1151; 1209; 1227
Selimov Renat	987; 1010	Shvartsburg Alexandre	199
Selzner Markus	698	Siciliano Pietro	366
Sendatzki Ann-Katrin	594	Sickmann Albert	436
Sendler Matthias	1039	Siebert Janusz	720
Sendyk Sandra	602	Siegers Wolter	1050

Sievert Julie	658	Smith Jason	772
Signorelli Paola	1138	Smith Richard D.	196
Sikorski Pawel	245	Smith Rob	552
Silverdale Monty	745	Snovida Sergei	438
Simair Monique	1058	Soares Renata	193
Simmonds MSJ	1041	Sobolev Anatoly P.	1105
Simmons Doug	614	Sobott Frank	12
Simon Aude	299	Söderström Martin	817
Simon Daniel	740	Sogawa Kazuyuki	276; 277
Simon Marc	307	Sojka Martha	1082
Simon Stéphanie	256; 440	Sojo L.	429
Simonato Manuela	39; 520; 744	Sokol Elena	1119
Sinclair Eleanor	745	Sokolová Jitja	697
Sinclair Ian	492	Solari Roberto	190
Singer Florian	454	Soleimani Ilnaz	377
Singer Heinz	184; 1049	Solouki Touradj	86
Singh Dilip Kumar	603	Solovyeva Elizaveta	223; 280
Singh Saranjit	603	Soltwisch Jens	87; 498; 632; 896
Singh Varoon	53	Sommer Ulf	1141
Sinisgalli Chiara	238	Somsen Govert W.	57; 529
Sinning Irmgard	116	Son Jin Gyeong	911
Sinz Andrea	118; 787	Sonego Roberto	1106
Siracusa Melania	1108	Song Jihyeon	1197
Siri Didier	1161	Song Ji-Inn	75
Široká Jitka	544	Song Lijiang	546
Siu Andy C. K.	938	Song Sanghoon	275
Siu Chi-Kit	447	Song Shuling	999
Siu K. W. Michael	129	Song Yoon Chae	1084
Siviero Andrea	1147	Sonnette Alexandre	350
Sivonen Kaarina	765	Soo Jung Hu	1110
Sjerps Rosa	1094	Sophie Siems	412
Sjöberg Per J. R.	1077	Sorokin Anatoly	40; 1151; 1209; 1227
Skalicka-Woźniak Krystyna	414; 525	Sosnowski Tori	291
Skehel Mark	781	Sottani Cristina	66
Skjødt Karsten	90	Southwick Katie	788
Sklorz Martin	82; 457	Souza Amanda	730; 731; 746; 961
Skraskova Karolina	772	Sovová Kristýna	182; 482
Skriba Anton	379; 1236	Spacil Zdenek	38; 432
Skrzypkowska Maria	720	Spadafora Natasha D.	875
Skudas Romas	517	Španěl Patrik	182; 482; 483; 944
Škvára Pavel	545; 605	Spasskaya Daria	405
Skwara O.	1045	Specht Harrison	69
Slavata Lukas	10	Spencer Jeremy P E	615
Slavov Nikolai	69	Spencer Simon	222
Sleno Lekha	141	Spengler Bernhard	259; 922
Slutsky Smith Elana	1140	Speranza Maurizio	448
Smietanski Mirosław	245	Sperling Michael	897; 914
Smit Elize	600	Spesyvyi Anatolii	182; 482; 483
Smith Andrew James	253; 271; 920; 921; 923	Speybrouck David	538
Smith Clive A.	331; 844	Spezia Riccardo	48; 112; 128
Smith David	482; 963	Spicka Jan	385
Smith Donald F.	78; 201	Spjuth Ola	1211

Spraggins Jeffrey	122	Strnad Miroslav	544; 733; 1213
Springorum Patrizia	787	Strojnik Lidija	853; 862
Squillaci Donato	846	Strupat Kerstin	259; 922
Srivastava Deepika	556	Stryjak Iga	698
Stacchini Paolo	613	Strzelecka Dominika	245
Staccioli Serena	837	Stuchfield Dale	173
Stacey Gary	81	Študent Vladimír	1126
Stafford George	1182	Studio Emanuel	66
Stahl-Zeng Jianru	979; 994; 1000; 1107	Stulik Jiri	425
Staines Stephen	35	Su Yong-yang	1085
Stanford Phoebe	695	Suarez-Fernandez Amanda	422
Starodubtseva Natalia	289	Subra Gilles	132
Starr Nichola	628	Suckau Detlev	559
Stashenko Elena	1043	Sudakov Mikhail	157
Štaud František	544	Sugiura Yuki	747; 919
Stavrianidi Andrey	819; 835	Sugiyama Eiji	747
Steacker Kelly	408	Sugiyama Masaaki	802
Stead Sara	198; 861	Sugiyama Masuyuki	494
Stecconi Tommaso	991; 1001; 1002; 1003	Suh Jung Ki	1109
Steckel Arnold	340; 650	Sujin Lee	639; 643
Stegmann Olmedillas Juan Luis	550	Sukau Detlev	1158
Steenvoorden Evelyne	904	Sukhikh Gennady	459
Stefanowicz Piotr	165; 1045; 1201; 1202	Sulek Karolina	767
Stefański Tomasz	522	Sulzer Philipp	478; 1080
Steffan Agostino	364; 703	Suman Michele	851; 1007
Steffensen Stine Krogh	388	Sun Xiuxuan	436
Steil Daniel	896	Sun Zhumei	157
Stein Stephen E.	89; 95; 167; 555; 757; 858; 1232	Supuran Claudiu T.	841; 846
Steinbeck Christoph	1211	Sushkov Oleg	736
Steinbrueck Clara	401	Sut Stefania	236
Steinhuber Bernd	935	Svatoš Aleš	434; 487; 626; 907
Stella Martina	253; 271; 921; 923	Svehla Elisabeth	1024
Stella Paola	984	Svirikova Anastasiya	913
Stellin Giovanni	744	Swainston Neil	1186
Steven Rory T.	147; 924; 1222	Syka John	956
Stevens Katherine	121	Sysoev Alexey	1179
Stiborek Marek	505	Szabo Istvan	308
Stienstra Rinke	278	Szakacs Gergely	740
Stöberl Anna-Sophie	729	Sze Andre	938
Stockett Mark	113; 309	Szesny Matthias	748; 1142
Stopka Sylwia A.	81	Szewczuk Zbigniew	165
Storti Matteo	520	Szymgel Łukasz	380
Stout Tom	563	Szpunar Joanna	501
Stramenga Arianna	991; 1001; 1002; 1003	Szuchman-Sapir Andrea	1140
Štrancar Ales	343	Szykula Katarzyna	773
Stratton Tim	823; 870	Szynkowska Małgorzata Iwona	664; 668; 673; 888
Streibel Thorsten	82; 181		
Streli Christina	913	T	
Strettoi Enrica	374	Tabet Jean-Claude	49
Strittmatter Nicole	898	Tack Cees	278
Strmeň Timotej	354	Tada Asa	636

Tahir Muhammad	436	Tchekhovskoi Dmitrii	89; 95; 1232
Tait Sabrina	1053; 1117	Tedesco Daniele	562
Takada Hidenori	623	Tejada-Casado Carmen	1165
Takada Yasuaki	494	Telensky Petr	38
Takahashi Eiichi	809	Tellström Verena	427
Takahashi Hidenori	151; 949	Temporini Caterina	541
Takahashi Hisashi	957	Ten Hacken Nick H. T.	265
Takano Eriko	1186	Tengattini Sara	541
Takano Shigetsugu	276	Tennagels Norbert	549
Takáts Zoltán	147; 193; 474; 740; 875; 898; 908; 924	Tenzer Stefan	412
Takayama Mitsuo	101	Teodori Elisabetta	840
Takei Chikako	369; 507; 508; 672; 891; 1192	Teppo Jaakko	444
Takeo Emi	919	Ter Laak Thomas	1050
Takeuchi Atsuko	389	Terada Takeshi	693
Takeuchi Aya	54; 330	Tereshkova Alesya	280
Talebpour Zahra	355	Terui Yasushi	609
Tamimi Sara	571; 1099	Tesler Larry	102
Tampé Robert	799	Tessaro Elias Paulo	827
Tan Gina	746	Tettamanti Silvia	271; 916
Tan Minjie	563	Thader Andreas	1020; 1023
Tan Nadia Suray	394	Thais Lino	884
Tan Warren K. Y.	749	Thakeow Prodpran	362
Tan Wei	11	Theisen Alina	305
Tanaka Koichi	151; 949; 957	Thelemann Jonas	74
Tanaka Koji	936	Thevis Mario	284; 349; 1183
Tanaka Toshiya	639	Thibault Pierre	44
Tang Wai Kit	447	Thiel Marcel	257
Tang Xiaofan	793	Thiel Sira	729
Tang Yang	683	Thissen Roland	950
Tang Zhixu	620; 621; 622; 1009	Thøgersen Ida B.	654
Taniguchi Moyu	727	Thomas Andreas	284; 1183
Tans Roel	278	Thomas Dennis G.	196
Taoka Masato	9	Thomas Dominique	1143
Tarasova Irina	223; 280; 405	Thomas Marion	750
Tart Jonathan	556	Thomas Mary	183
Tascon Marcos	53; 65; 383; 528	Thomas Sebastien	538; 1087
Tassinari Roberta	1117	Thomas Spencer A.	147; 924; 1222
Tata Alessandra	191	Thompson Christopher	80
Tatam Ralph P.	35	Thon Vojtech	432
Tatituri R	1150	Thorn James	209
Taubert Anja	922	Thorn Jim	312; 529; 647; 698
Tauler Ferre Romà	1136	Thorsteinsdóttir Margrét	396
Tautenhahn Ralf	730; 731; 746; 823; 870	Thorsteinsdóttir Unnur Arna	396
Tavani Fabiana	824; 838	Tian Xiaobo	794
Taverna Domenico	443	Tian Yuan	114
Tavoloni Tamara	991; 1001; 1002; 1003; 1234	Tiberi Paolo	925
Taxer Thomas	104; 304	Tielens Alexander	109
Taylor Adam	147; 804; 924; 1222	Tierny Dominique	458
Taylor Adrian	660	Tietz Drew	1138
Taylor Michael	191; 490	Tikhonov Georgy	311; 575
Taylor Phil	595	Tilemann Lisa	518
		Timoumi Laila	1157

Tintaru Aura 361; 547
 Tinti Michele 437
 Tisato Francesco 105
 Toffoli Giuseppe 513
 Togna Anna Rita 671
 Tokarski Caroline 225
 Tölgyesi László 604
 Tolpina Miriam 591
 Tomaniova Monika 598
 Tomati Valeria 403
 Tomonaga Takeshi 263; 510; 636
 Tonazzini Ilaria 423
 Tonge Michael 556
 Tonidandel Loris 857; 1011
 Torma Frank A. 445
 Toro Botros 187
 Torres Sonia 751
 Torta Federico 1148
 Tortorella Sara 925
 Tóth Gábor 1028
 Touboul David 1044; 1144
 Tovar Trenton 473
 Town James 509
 Trainotti Debora 1011
 Traldi Pietro 571; 1037; 1099; 1242
 Tamarin Anna 281
 Tran Bao 883
 Tranvik Lars J. 1077
 Tratkiewicz Ewa 830
 Traub Heike 206
 Travelli Cristina 512
 Treit Peter 286
 Tremintin Guillaume 559; 1158
 Tremoli Elena 753
 Trengove Robert 719
 Treu Axel 910
 Treutler Hendrik 709; 713
 Trevisiol Stéphane 895
 Trevitt Adam 21; 125
 Triolo Antonio 133; 824; 831; 837; 838
 Trisna Meyci 97
 Trivedi Drupad 745
 Troensegaard Nielsen Kim 903
 Troiani Anna 946; 947; 948; 951
 Troiani Sonia 548
 Trojsi Giorgio 1244
 Tsai Shang-Ting 1021
 Tsai Ya-Ting 1113
 Tseng Mei-Chun 1029
 Tsipi Despina 1109
 Tsou Hsin-Pai 1133
 Tsuchida Masako 1201

Tsugoshi Takahisa 353
 Tsuji Makoto 753
 Tsybin Oleg Y. 962
 Tsybin Yury O. 159; 631; 958; 962
 Tung M. 429
 Tupinier Jérôme 350
 Turan Yucel Nazli 821
 Turco Francesca 863
 Turecek Frantisek 20
 Turiák Lilla 686; 1028
 Turk John 1150
 Turkina Maria 258; 430
 Turkseven Seyma 193
 Turner Claire 367; 963
 Turner Matthew 773
 Turnu Linda 753
 Turyanskaya Anna 913

U

Ube Yuko 623
 Ubhi Baljit K. . . 708; 723; 724; 726; 741; 754; 1120
 Ubukata Masaaki 1228
 Uchimura Tatsuo 623
 Ueda Hiroki 639; 643
 Ueda Yoshihisa 1228
 Uetrecht Charlotte 800; 805
 Uhl Jenny 80
 Uhrová Martina 989
 Ujma Jakub 173
 Ulisses Tega David 317; 689
 Ulyanenkov Alex 1220
 Ul'yanovskii Nikolay 1064
 Umekage So 693
 Undas Anna 211
 Unger Tamar 798
 Uppal Annu 521; 535; 785; 1173
 Upton Rosie 67
 Uray Katalin 340; 650
 Urbánek Lubor 544
 Urteil Jonathan 1020
 Usai Donatella 1242
 Usen Ime C 356

V

Vadi Alessandro 796
 Vaglio Sara 372
 Vaikkinen Anu 444; 770
 Vain Thomas 733
 Vaitsekhovich Viktoryia 550
 Valera Stefano 1048
 Valis Karel 10
 Valledor González Luis 419

Vallet Marine	907	Ventura Rosa	554
Valsasina Barbara	548	Venturelli Matteo	1018
Valsesia Armand	250	Verbeke Lynn	426; 431
van Agthoven Maria	158; 964; 1230	Verbick Newton	189
van Bommel Maarten	226	Verdegaal Suzan	97
van Dam Michael	160	Vereyken Liesbeth	131
van de Bram Put	57	Verhaert Peter Dem	269
van de Merbel N.C.	265	Verma Alka	994
van de Waterbeemd Michiel	156	Vernhes Emeline	117
van den Berg Klaas Jan	226	Vertes Akos	81
van der Lee Martijn	211	Verzeaux Laurie	1041
van der Linde Christian	20; 104; 107; 301; 304	Vetere Alessandro	1229
van der Meer Bram	586; 1126	Veyel Daniel	1172
Van der Meer Jeroen A.	586	Viani Mario	497
Van der Rest Guillaume	1163	Vidale Massimo	230
van der Schans Marcel J.	586	Vidova Veronika	432
van der Schee Marc	552	Viehoff Maria	897
Van Duijneveldt Jeroen	1125	Viel Diane	407
van Geenen Fred	127	Viet Laura	876
van Gool Alain	278	Vigne Sébastien	350
Van Hoeck Els	163; 1171	Vikström Elena	430
van Outersterp Rianne	725	Villar Garea Ana	549
van Pijkeren Alienke	650	Villette Claire	890
van Rooij Arno	278	Vimer Shay	115; 777; 798
van Ruth Saskia	135	Vinati Simone	247
van Tol Helena	563	Vincendet Jean-Baptiste	754
Van Veelen Peter	921	Vincenti Flaminia	669; 670; 671; 1098
van Wasen Sebastian	484	Vincenti Marco	372
VanAernum Zachary	778	Viner Rosa	433; 781; 792; 799
Vane Christopher	183	Viola Andrea	316
Vanni Ester	16	Vissers Hans	1172
Vanninen Paula	143; 669; 817	Vissers Johannes	1190
Vannucchi Francesca	1088	Vitek Olga	17; 426
Vardanega Renata	1036	Vittori Elisa	1102
Vargas Pérez Marta	979	Vittori Sauro	24; 1096; 1102; 1206
Vargas Thais	996	Vivaldi Barbara	970; 972
Varmuza Kurt	691	Vladimirov Gleb	157; 928
Vasiljevic Tijana	383; 485; 486	Voeten Robert	57
Vasilopoulou Catherine	767	Vogel Heiko	434
Vasilyeva Alexandra	651	Vogel Ulla	248
Vázquez de Ágredos Pascual Maria Luisa	1242	Voinov Valery	676; 1182
Vedder Sven	1101	Vojs Marian	605
Vedovelli Luca	39; 520; 744	Vojs Staňová Andrea	453; 545; 605; 1089
Veggi Daniele	1034	Vojta Simon	587
Veglia Fabrizio	753	Völklein Friedemann	28; 536
Vékey Károly	686; 1028	Volland Hervé	440
Veličković Dusan	81	Volmer Dietrich A.	456; 484; 1181
Velickovic Suzana	839	Volny Michael	385
Veljkovic Filip	839	Voloaca Oana	929
Venier Marta	616	von Gunten Urs	601
Ventura Emiliano	144	von Woedtke Thomas	636; 653
Ventura Giovanni	316	Vonka Petr	728

Vorm Ole 286
 Vorng Jean-Luc 147
 Vorobyev Alexander 40; 1227
 Vosse Christian 1129
 Vrána David 1126
 Vreeken Rob J. 63; 131; 584; 629; 832
 Vrhovsek Urska 705
 Vrkoslav Vladimír 354
 Vu Huu Khanh 799
 Vughs Dennis 1050
 Vukelić Željka 1030; 1178
 Vydelingum Seeven 134

W

Wada Hiroshi 564
 Wada Motoi 949
 Wagner Roger 996
 Wåhlén Karin 283
 Wahlsten Matti 765
 Waidelich Dietmar 413
 Wakelam Michael 190
 Waliczek Mateusz 165
 Wallace Federico 346
 Wallace William 95; 1221
 Walpurgis Katja 284
 Walter Kerstin 910
 Walter Ondine 250
 Walter Stefan 404
 Walton Chris 35
 Wang Alexandre 660
 Wang Bing 551
 Wang Huixin 221
 Wang Kelin 86
 Wang QianQian 579
 Wang Tianfang 389
 Wang Wei 359; 1085
 Wang Wenliang 357; 607
 Wang Wenning 567
 Wang Yan 52
 Wang Yuning 624
 Ware Rebecca L. 78
 Warminski Marcin 245
 Warren John 445
 Warskulat Anne-Christin 434
 Warszawsky Shira 777
 Warzok Ulrike 106
 Washiyama Yuya 328
 Wasilewska Alicja 572
 Wąsowicz Marcin 383
 Watanabe Hiroshi 715
 Watanabe Jun 381
 Watolla Carsten 634
 Watts Andrew 539
 Watts Peter 887
 Webb Ian K 196
 Webb Sally 755; 1145
 Weber Darren 285
 Weber Ronja 454
 Wei Guanyi 357
 Wei Juan 683
 Weis Patrick 296
 Weisbrod Chad R. 78
 Weiss Victor U. 800
 Weissflog Jerrit 487
 Welander Jenny 430
 Welkie David 820
 Welte Stefan 549
 Weltmann Klaus-Dieter 636; 653
 Wende Kristian 636; 653
 Wendt Jürgen 203
 Wendt Karin 854
 Wendt Michael 106
 Wenk Markus 1148
 Wenske Sebastian 636; 653
 Wenzel Gabi 299
 Werner Hauke 412
 Wesdemiotis Chrys 56
 Wessels Hans 278
 Wessjohann Ludger 594
 West Andy 70
 West Brandi 153; 690
 Westphall Michael 1184
 Wetzel Dorothee 618; 656
 Wevers Ron 294; 725
 White Robert L. 684
 Whitelegge Julian 608
 Whitmore Christopher 358
 Widdig Anja 337
 Wieczorek Maciej 830
 Wiederhold Katrin 781
 Wiegand Pascal 28; 281; 536
 Wielogorska Ewa 852
 Wielsch Natalie 434
 Wierzbicka Magdalena 165; 1045
 Wijyantha-Kahagala-Gamage Upul 356
 Wilharm Thomas 203
 Wilkes Heinz 700
 Wilkins Charles 488
 Willets Matt 649
 Willetts Matthew 1230
 Willms J. Alexander 816
 Wilmott James S. 267
 Wilson Derek 172; 789
 Wilson Ian 719

Wilson John 446
 Wilson Richard 772
 Wingfield Jon 492
 Wingfield Jonathan 556
 Winkler Klaus 478
 Winkler Robert 96; 489; 630
 Winnie Hung Y.-L. 85
 Winter Dominic 397; 421; 640
 Wisthaler Armin 1080
 Witt Matthias 183; 204; 549; 890
 Wittig Sabine 795; 801
 Wojciechowski Kamil 1025
 Wojcik Roza 196
 Wojnicz Aneta 526; 550
 Wolczanski Grzegorz 1202
 Wolf Jana 781
 Wolf Michael 113
 Wolkoff Peder 1062
 Wollscheid Bernd 436
 Wolrab Denise 1126
 Wolski Witold E. 562
 Wood Troy 742
 Woodall Daniel 195; 1035; 1175
 Wooding Madelien 29
 Woods Lucy 275
 Woolman Michael 191; 490
 Wootton A. Christopher 1230
 Wootton Christopher A. 75; 91; 158; 959; 964
 Worbs Sylvia 440
 Wouters Elroy 1184
 Wrona Mark 573; 1217
 Wu Jingcun 337; 992; 1060; 1091
 Wu Megan 191; 490
 Wu Qinghao 114
 Wu Qiong 959
 Wu Ri 1185
 Wu Tong 407
 Wu Wei 563
 Wuhler Manfred 166
 Wyeth Russell 389
 Wysocki Vicki 75; 778

X

Xenopoulos Alex 517
 Xi C. 429
 Xia Yu 126
 Xie Fang 925
 Xie Yongming 337
 Xiong Lei 1120
 Xu Fuxing 323
 Xu Jiang 359
 Xu Shuting 450; 965

Xu Yun 205; 745
 Xuan Yue 410; 435; 436

Y

Yablokov Artur 928
 Yadav Avinash 679
 Yajima Daisuke 891
 Yalcin Talat 952
 Yamada Masaki 1122
 Yamada Naohito 756
 Yamaguchi Shinichi 930
 Yamamoto Atsushi 1004
 Yamamoto Kazuki 639; 643
 Yamanaka Sakino 276
 Yamauchi Yoshio 9
 Yamazaki Makoto 930
 Yamazaki Yuji 623
 Yamazaki Yuzo 802; 930
 Yan Bin 147; 804; 924; 1222
 Yan Cunyu 1186
 Yan Hao 805
 Yan Jianfei 1060
 Yan Rongxin 1090
 Yan Xinjian 89
 Yanes Oscar 905
 Yang Dorothy 604
 Yang Mo 464; 466
 Yang Pengyi 399
 Yang Pengyuan 1231
 Yang Xiaoyu 757
 Yang Xiaoyue 126
 Yang Yi 1231
 Yang Yongliang 1071
 Yang Zhibo 68
 Yano Maho 277
 Yao Zhong-Ping 447
 Yasueda Hisashi 693
 Yates John R. 164
 Ye Josh 820; 967
 Ye Sha Josh 1091
 Ye Xiwen 620; 621; 622; 1009
 Yefchak George 619
 Yeo Woon-Seok 1195; 1197
 Yi Bongyoon 332; 464; 466
 Yi Long 11
 Yi Sung 1066
 Yildirim Keleş Z. Zeynep 511
 Yim Yong-Hyeon 1109
 Yin Lois 134
 Yli-Kauhaluoma Jari 814
 Yokoi Yasuto 1145; 1156
 Yokoyama Ryo 785

Yona Meital 798
 Yong-Woon Shin 1110
 Yoo Jong Shin 1217
 Yoon Dahye 766
 Yoon Sohee 911
 Yoon Sung Hwan 1123
 Yoon Young-Ran 714; 718
 Yoshinari Kiyomi 609
 Yoshitomi Hideyuki 276
 Yoshiya Kazumi 1216
 Yoshizawa Kenichi . 369; 507; 508; 672; 891; 1192
 Young Phoebe 241; 685
 Young Reuben 192
 Yu Xiangying 447
 Yu Ying Qing 517; 574
 Yuan Gu 11
 Yuan Jia 149
 Yuan Xianglong 357; 359; 1215
 Yue Siyao 1167
 Yuneva Mariia O. 147; 924
 Yunin Maxim 892
 Yurina Lyubov' 651
 Yuste Córdoba Francisco Javier 1095

Z

Zaar Oscar 934
 Zabka Jan 110
 Zabrouskov Vlad 269; 428; 731; 823
 Zaitso Kei 218
 Zakharova Natalia 292
 Załęcka Karolina 668
 Zambrano Jesus 1043
 Zamfir Alina Diana 1030; 1178
 Zamith Sebastien 688
 Zamora Ismeal 1217
 Zampieri Davila 1111
 Zanella Veronica 1093
 Zang JieChao 677
 Zangl Rene 790; 1169; 1187
 Zangrando Roberta 1093
 Zapal Jakub 1236
 Zara Vincenzo 366
 Zarrine-Afsar Arash 191; 490
 Zarrineh Mahshid 377
 Zazza Costantino 448
 Zeegers Guido Paul 120
 Zelei Georgina 470
 Zeller Martin 687
 Zeman Florence 885
 Zenaidee Muhammad 154
 Zenkevich Igor 1072
 Zenobi Renato 51; 74; 120; 216;
 454; 475; 520; 729; 811
 Zerega Yves 463
 Zetterberg Henrik 249; 954
 Zettergren Henning 113; 309
 Zhai Lihua 357; 359; 607; 1215
 Zhang Chao 932
 Zhang Fan 1135
 Zhang Feiyu 171
 Zhang Hong 134
 Zhang Li 620
 Zhang Li 621
 Zhang Lihua 1026
 Zhang Lulu 11
 Zhang Pingyu 91
 Zhang Qifeng 190
 Zhang QingMin 1215
 Zhang Rutan 758
 Zhang Tao 794
 Zhang Terry 788
 Zhang Yiding 25
 Zhang Zheng 1232
 Zhao Fei 1009
 Zhao Shuang 748
 Zhao Wenbo 1093
 Zhen Zhenpeng 967
 Zheng Yufang 1046
 Zherebker Alexander 178; 468; 928
 Zhokhov Sergey 591
 Zhong Feng 614
 Zhou Hu 436
 Zhou Jiang 11
 Zhou Mowei 658
 Zhou Xiangdong 337
 Zhou Yue 435; 436
 Zhu Hongwen 436
 Zhu Mingli 967
 Zhu Zhentao 11
 Zhvansky Evgeny 40; 1151; 1209; 1227
 Zibara Kazem 148
 Zifferero Andrea 230
 Zigon Dusan 527
 Ziller Luca 862
 Zimens Marina 818
 Zimmerli Christian 17
 Zimmermann Ralf 82; 181; 203; 457
 Zimmermann Stefan 943; 1157; 1176
 Zimnicka Magdalena 449; 572
 Ziolkowski Robert 524
 Ziskind Michael 458
 Zocolo Guilherme Julião 1111

Zompanti Alessandro	670	Žukauskaitė Asta	733
Zoratto Samuele	800	Žukowicz Wioletta	542
Zsuga Miklós	59; 470; 477; 503	Zuo Xuan	893
Zubarev Roman A.	398; 682	Zuo Zhicheng	567
Zuberovic Muratovic Aida	1012	Zurowietz Annika	901
Zuidhof Martin	233	Zvereva Irina	892
Zuilhof Han	127	Zymak Illia	110



ISBN 9788890738852